

multichannel* systems

McsUsbNet.dll

Version 5.0.11

Multi Channel Systems MCS GmbH
Aspenhastrasse 21
72770 Reutlingen
Germany
Fon +49-71 21-90 92 5 - 0
Fax +49-71 21-90 92 5 -11
info@multichannelsystems.com
www.multichannelsystems.com

Generated by Doxygen 1.8.18

1 McsUsbNet.dll for MCS USB devices	1
1.1 Introduction	1
1.2 System requirements	1
2 Other MCS USB devices	2
2.1 MCS-USB-Sw2to64 device	2
2.2 The MCS FluidControl Device	2
2.2.1 Introduction	2
2.2.2 Access to the FluidControl device	3
3 MC_Card and MCS-USB-ME(A) Data Acquisition Devices (DAQ)	3
4 The MCS Robo Device	4
4.1 Introduction	4
5 STG200x & STG400x Series Stimulus Generators	4
5.1 Introduction	4
5.2 Download mode	5
5.2.1 Memory Layout and Trigger Setup	5
5.3 Streaming mode	7
5.3.1 Memory Layout and Trigger Setup	8
6 Namespace Index	10
6.1 Namespace List	10
7 Hierarchical Index	10
7.1 Class Hierarchy	10
8 Class Index	15
8.1 Class List	15
9 Namespace Documentation	21
9.1 Mcs Namespace Reference	21
9.2 Mcs::Usb Namespace Reference	21
9.2.1 Enumeration Type Documentation	25
9.2.2 Function Documentation	25
10 Class Documentation	27
10.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference	27
10.1.1 Member Data Documentation	27
10.2 BatteryState Class Reference	28
10.2.1 Property Documentation	28
10.3 BesselFilterHighPassNet Class Reference	28
10.3.1 Constructor & Destructor Documentation	29
10.4 BesselFilterLowPassNet Class Reference	29
10.4.1 Constructor & Destructor Documentation	29

10.5 ButterworthFilterHighPassNet Class Reference	29
10.5.1 Constructor & Destructor Documentation	30
10.6 ButterworthFilterLowPassNet Class Reference	30
10.6.1 Constructor & Destructor Documentation	30
10.7 CChannelTestDeviceNet Class Reference	31
10.7.1 Constructor & Destructor Documentation	31
10.7.2 Member Function Documentation	31
10.8 CCMOSMea_FunctionNet Class Reference	32
10.8.1 Constructor & Destructor Documentation	34
10.8.2 Member Function Documentation	34
10.9 CCMOSMeaDeviceNet Class Reference	43
10.9.1 Constructor & Destructor Documentation	44
10.9.2 Member Function Documentation	44
10.9.3 Property Documentation	45
10.10 CCreateFilterNet Class Reference	46
10.10.1 Constructor & Destructor Documentation	46
10.10.2 Member Function Documentation	47
10.10.3 Property Documentation	47
10.11 CDacCalibrationFunctionNet Class Reference	48
10.11.1 Detailed Description	48
10.11.2 Constructor & Destructor Documentation	48
10.11.3 Member Function Documentation	49
10.12 CDacqGroupChannelGenericSelectionNet Class Reference	49
10.12.1 Constructor & Destructor Documentation	50
10.13 CDacqGroupChannelSelectionNet Class Reference	50
10.13.1 Constructor & Destructor Documentation	50
10.14 CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, Dacq↔ GroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference	51
10.14.1 Constructor & Destructor Documentation	51
10.14.2 Member Function Documentation	51
10.15 CDeviceGroupChannelInfoGenericNet Class Reference	53
10.15.1 Constructor & Destructor Documentation	54
10.16 CDeviceGroupChannelInfoNet Class Reference	54
10.16.1 Constructor & Destructor Documentation	54
10.17 CDeviceGroupChannelInfoSCUNet Class Reference	54
10.17.1 Constructor & Destructor Documentation	55
10.18 CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet > Class Tem- plate Reference	55
10.18.1 Constructor & Destructor Documentation	55
10.18.2 Member Data Documentation	55
10.19 CDeviceGroupChannelInfoW2100Net Class Reference	56
10.19.1 Constructor & Destructor Documentation	56
10.20 CDigOutStimulatorFunctionNet Class Reference	56

10.20.1 Detailed Description	57
10.20.2 Constructor & Destructor Documentation	57
10.20.3 Member Function Documentation	58
10.21 CEncapsulatorDeviceNet Class Reference	60
10.21.1 Detailed Description	61
10.21.2 Constructor & Destructor Documentation	61
10.21.3 Member Function Documentation	61
10.22 CExternDTesterDeviceNet Class Reference	61
10.22.1 Detailed Description	62
10.22.2 Constructor & Destructor Documentation	62
10.22.3 Member Function Documentation	62
10.23 CFilterCoefficientsNet Class Reference	63
10.23.1 Constructor & Destructor Documentation	64
10.23.2 Member Function Documentation	65
10.23.3 Property Documentation	65
10.24 CFilterConfigurationNet Class Reference	67
10.24.1 Constructor & Destructor Documentation	67
10.24.2 Member Function Documentation	67
10.25 CFilterConfigurationRegisterNet Class Reference	68
10.25.1 Constructor & Destructor Documentation	69
10.25.2 Member Function Documentation	69
10.26 CFilterPropertyNet Class Reference	70
10.26.1 Constructor & Destructor Documentation	70
10.26.2 Member Function Documentation	71
10.26.3 Property Documentation	71
10.27 CFluidControlDeviceNet Class Reference	72
10.27.1 Detailed Description	73
10.27.2 Constructor & Destructor Documentation	73
10.27.3 Member Function Documentation	73
10.27.4 Property Documentation	78
10.28 CFYIDeviceNet Class Reference	78
10.28.1 Detailed Description	78
10.28.2 Constructor & Destructor Documentation	78
10.28.3 Property Documentation	79
10.29 CGenericDevelopDeviceNet Class Reference	79
10.29.1 Detailed Description	86
10.29.2 Constructor & Destructor Documentation	86
10.29.3 Member Function Documentation	86
10.30 CGilsonDeviceNet Class Reference	96
10.30.1 Detailed Description	97
10.30.2 Constructor & Destructor Documentation	97
10.30.3 Member Function Documentation	97

10.30.4 Member Data Documentation	98
10.31 CHiClampDeviceNet Class Reference	98
10.31.1 Detailed Description	99
10.31.2 Constructor & Destructor Documentation	99
10.31.3 Property Documentation	99
10.32 CHLADacqNet Class Reference	99
10.32.1 Constructor & Destructor Documentation	100
10.33 CHLADeviceNet Class Reference	100
10.33.1 Detailed Description	100
10.33.2 Constructor & Destructor Documentation	100
10.33.3 Property Documentation	101
10.34 CMcsUsbDacqNet::CHWInfo Class Reference	101
10.34.1 Detailed Description	101
10.34.2 Constructor & Destructor Documentation	101
10.34.3 Member Function Documentation	102
10.35 CIntanMea_FunctionNet Class Reference	103
10.35.1 Constructor & Destructor Documentation	104
10.35.2 Member Function Documentation	104
10.36 CInterfaceboardFunctionNet Class Reference	105
10.36.1 Detailed Description	106
10.36.2 Constructor & Destructor Documentation	106
10.36.3 Member Function Documentation	107
10.37 CLIH3DeviceNet Class Reference	107
10.37.1 Detailed Description	109
10.37.2 Constructor & Destructor Documentation	109
10.37.3 Member Function Documentation	109
10.37.4 Property Documentation	114
10.38 CMcsBus_AxisParametersNet Class Reference	114
10.38.1 Constructor & Destructor Documentation	115
10.38.2 Member Function Documentation	115
10.39 CMcsBus_ExtensionNet Class Reference	116
10.39.1 Constructor & Destructor Documentation	116
10.39.2 Member Function Documentation	117
10.40 CMcsBus_FYIExtensionNet Class Reference	117
10.40.1 Constructor & Destructor Documentation	117
10.40.2 Member Function Documentation	118
10.41 CMcsBus_MotorControlNet Class Reference	119
10.41.1 Constructor & Destructor Documentation	122
10.41.2 Member Function Documentation	122
10.42 CMcsBus_SensorNet Class Reference	136
10.42.1 Constructor & Destructor Documentation	138
10.42.2 Member Function Documentation	138

10.43 CMcsBus_TempSensorNet Class Reference	146
10.43.1 Constructor & Destructor Documentation	147
10.43.2 Member Function Documentation	147
10.44 CMcsBus_VoltageModeNet Class Reference	148
10.44.1 Constructor & Destructor Documentation	149
10.44.2 Member Function Documentation	150
10.45 CMcsBusNet Class Reference	153
10.45.1 Constructor & Destructor Documentation	153
10.45.2 Member Function Documentation	154
10.46 CMcsUsbDacqNet Class Reference	157
10.46.1 Detailed Description	162
10.46.2 Constructor & Destructor Documentation	163
10.46.3 Member Function Documentation	163
10.46.4 Member Data Documentation	204
10.46.5 Property Documentation	204
10.46.6 Event Documentation	205
10.47 CMcsUsbDeviceStatePushFunctionNet Class Reference	205
10.47.1 Constructor & Destructor Documentation	205
10.47.2 Member Function Documentation	206
10.47.3 Event Documentation	206
10.48 CMcsUsbDeviceStatePushNet Class Reference	206
10.48.1 Constructor & Destructor Documentation	206
10.48.2 Member Function Documentation	207
10.48.3 Event Documentation	207
10.49 CMcsUsbFactoryNet Class Reference	207
10.49.1 Constructor & Destructor Documentation	209
10.49.2 Member Function Documentation	209
10.49.3 Member Data Documentation	215
10.50 CMcsUsbFunctionNet Class Reference	215
10.50.1 Constructor & Destructor Documentation	216
10.50.2 Member Function Documentation	216
10.50.3 Member Data Documentation	217
10.51 CMcsUsbFunctionPointerContainer Class Reference	217
10.52 CMcsUsbListEntryNet Class Reference	217
10.52.1 Detailed Description	218
10.52.2 Constructor & Destructor Documentation	218
10.52.3 Member Function Documentation	218
10.52.4 Property Documentation	220
10.53 CMcsUsbListNet Class Reference	221
10.53.1 Detailed Description	222
10.53.2 Constructor & Destructor Documentation	222
10.53.3 Member Function Documentation	222

10.53.4 Property Documentation	223
10.53.5 Event Documentation	223
10.54 CMcsUsbNet Class Reference	224
10.54.1 Detailed Description	228
10.54.2 Constructor & Destructor Documentation	228
10.54.3 Member Function Documentation	228
10.54.4 Member Data Documentation	240
10.54.5 Property Documentation	245
10.55 CMcsUsbPointerContainer Class Reference	245
10.56 CMEA2100x256FunctionNet Class Reference	245
10.56.1 Detailed Description	245
10.56.2 Constructor & Destructor Documentation	245
10.56.3 Member Function Documentation	246
10.57 CMeaAudioFunctionNet Class Reference	247
10.57.1 Constructor & Destructor Documentation	247
10.57.2 Member Function Documentation	248
10.58 CMeaCleanDeviceNet Class Reference	250
10.58.1 Detailed Description	251
10.58.2 Constructor & Destructor Documentation	251
10.58.3 Member Function Documentation	252
10.59 CMeaCoatDeviceNet Class Reference	254
10.59.1 Detailed Description	255
10.59.2 Constructor & Destructor Documentation	256
10.59.3 Member Function Documentation	256
10.60 CMeaDeviceNet Class Reference	260
10.60.1 Detailed Description	261
10.60.2 Constructor & Destructor Documentation	261
10.60.3 Member Function Documentation	262
10.60.4 Property Documentation	267
10.61 CMeaDigitalDataFunctionNet Class Reference	268
10.61.1 Constructor & Destructor Documentation	268
10.61.2 Member Function Documentation	269
10.62 CMeaFeedbackFunctionNet Class Reference	270
10.62.1 Constructor & Destructor Documentation	271
10.62.2 Member Function Documentation	271
10.63 CMeaImpedanceDeviceNet Class Reference	274
10.63.1 Constructor & Destructor Documentation	274
10.63.2 Member Function Documentation	275
10.64 CMeasureTableDeviceNet Class Reference	275
10.64.1 Detailed Description	276
10.64.2 Constructor & Destructor Documentation	276
10.64.3 Property Documentation	276

10.65 CMeaSwitchDeviceNet Class Reference	276
10.65.1 Detailed Description	277
10.65.2 Constructor & Destructor Documentation	277
10.65.3 Member Function Documentation	277
10.66 CMeaUSBDeviceNet Class Reference	278
10.66.1 Detailed Description	279
10.66.2 Constructor & Destructor Documentation	279
10.67 CMeFunctionNet Class Reference	279
10.67.1 Detailed Description	280
10.67.2 Constructor & Destructor Documentation	280
10.67.3 Member Function Documentation	280
10.68 CMultiBatteryChargerDeviceNet Class Reference	281
10.68.1 Detailed Description	282
10.68.2 Constructor & Destructor Documentation	282
10.68.3 Member Function Documentation	282
10.69 CMultiwellCallbackFunctionNet Class Reference	288
10.69.1 Detailed Description	289
10.69.2 Constructor & Destructor Documentation	289
10.69.3 Member Function Documentation	289
10.69.4 Event Documentation	291
10.70 CMultiwellDeviceNet Class Reference	291
10.70.1 Detailed Description	292
10.70.2 Constructor & Destructor Documentation	293
10.70.3 Member Function Documentation	293
10.71 CMultiwellOptoStimFunctionNet Class Reference	297
10.71.1 Detailed Description	297
10.71.2 Constructor & Destructor Documentation	297
10.71.3 Member Function Documentation	298
10.72 CNF_GenDeviceNet Class Reference	302
10.72.1 Constructor & Destructor Documentation	302
10.72.2 Member Function Documentation	302
10.73 COctoPotDeviceNet Class Reference	302
10.73.1 Constructor & Destructor Documentation	303
10.73.2 Member Function Documentation	303
10.74 COkuvisionStimulatorDeviceNet Class Reference	306
10.74.1 Constructor & Destructor Documentation	307
10.74.2 Member Function Documentation	307
10.75 CPatchServerDeviceNet Class Reference	310
10.75.1 Detailed Description	310
10.75.2 Constructor & Destructor Documentation	310
10.75.3 Property Documentation	310
10.76 CPathIdentDeviceNet Class Reference	311

10.76.1 Constructor & Destructor Documentation	311
10.76.2 Member Function Documentation	311
10.77 CPedoterDeviceNet Class Reference	312
10.77.1 Detailed Description	312
10.77.2 Constructor & Destructor Documentation	312
10.77.3 Member Function Documentation	312
10.78 CPeristalticPumpDeviceNet Class Reference	313
10.78.1 Detailed Description	313
10.78.2 Constructor & Destructor Documentation	314
10.78.3 Property Documentation	314
10.79 CPgaDeviceNet Class Reference	314
10.79.1 Constructor & Destructor Documentation	315
10.79.2 Member Function Documentation	315
10.80 CPositionIIDeviceNet Class Reference	316
10.80.1 Detailed Description	317
10.80.2 Constructor & Destructor Documentation	317
10.80.3 Member Function Documentation	317
10.80.4 Property Documentation	319
10.81 CPositionImpDeviceNet Class Reference	319
10.81.1 Detailed Description	320
10.81.2 Constructor & Destructor Documentation	320
10.81.3 Member Function Documentation	321
10.82 CPPCDeviceNet Class Reference	322
10.82.1 Constructor & Destructor Documentation	323
10.82.2 Property Documentation	323
10.83 CPPCFunctionNet Class Reference	324
10.83.1 Detailed Description	325
10.83.2 Constructor & Destructor Documentation	325
10.83.3 Member Function Documentation	325
10.84 CPPS_DeviceNet Class Reference	330
10.84.1 Constructor & Destructor Documentation	331
10.84.2 Property Documentation	331
10.85 CPPS_FunctionNet Class Reference	331
10.85.1 Constructor & Destructor Documentation	332
10.85.2 Member Function Documentation	333
10.86 CPPSDeviceNet Class Reference	336
10.86.1 Detailed Description	336
10.86.2 Constructor & Destructor Documentation	336
10.87 CProgramPressureCurveNet Class Reference	336
10.87.1 Detailed Description	337
10.87.2 Constructor & Destructor Documentation	337
10.87.3 Member Function Documentation	337

10.88 CPulseGeneratorFunctionNet Class Reference	338
10.88.1 Detailed Description	339
10.88.2 Constructor & Destructor Documentation	339
10.88.3 Member Function Documentation	339
10.89 CRadioControlledDevicesNet Class Reference	341
10.89.1 Constructor & Destructor Documentation	342
10.89.2 Member Function Documentation	342
10.90 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference	343
10.90.1 Constructor & Destructor Documentation	343
10.90.2 Member Function Documentation	343
10.90.3 Member Data Documentation	343
10.91 CRetinaLedDeviceNet Class Reference	344
10.91.1 Constructor & Destructor Documentation	344
10.91.2 Member Function Documentation	344
10.92 CRFFFunctionNet Class Reference	346
10.92.1 Detailed Description	346
10.92.2 Constructor & Destructor Documentation	347
10.92.3 Member Function Documentation	347
10.93 CRobo_FYIProgram_FunctionNet Class Reference	350
10.93.1 Constructor & Destructor Documentation	351
10.93.2 Member Function Documentation	351
10.94 CRobo_FYITemp_FunctionNet Class Reference	352
10.94.1 Constructor & Destructor Documentation	352
10.94.2 Member Function Documentation	353
10.95 CRoboDacqNet Class Reference	354
10.95.1 Constructor & Destructor Documentation	356
10.95.2 Member Function Documentation	356
10.96 CRoboDeviceNet Class Reference	367
10.96.1 Detailed Description	370
10.96.2 Constructor & Destructor Documentation	370
10.96.3 Member Function Documentation	371
10.96.4 Member Data Documentation	377
10.96.5 Property Documentation	381
10.96.6 Event Documentation	381
10.97 CRoboFluidDeviceNet Class Reference	381
10.97.1 Constructor & Destructor Documentation	382
10.97.2 Member Function Documentation	382
10.97.3 Member Data Documentation	384
10.97.4 Property Documentation	384
10.98 CRoboInjectDeviceNet Class Reference	385
10.98.1 Detailed Description	385
10.98.2 Constructor & Destructor Documentation	385

10.99 CRoboocyte2DeviceNet Class Reference	385
10.99.1 Detailed Description	386
10.99.2 Constructor & Destructor Documentation	386
10.99.3 Member Function Documentation	386
10.100 CRoboStatorDeviceNet Class Reference	387
10.100.1 Constructor & Destructor Documentation	388
10.100.2 Member Function Documentation	388
10.100.3 Property Documentation	393
10.101 CSafeISDeviceNet Class Reference	393
10.101.1 Detailed Description	394
10.101.2 Constructor & Destructor Documentation	394
10.101.3 Member Function Documentation	394
10.101.4 Property Documentation	395
10.102 CSCUDacqGroupChannelSelectionNet Class Reference	396
10.102.1 Constructor & Destructor Documentation	396
10.103 CSCUFunctionNet Class Reference	397
10.103.1 Detailed Description	399
10.103.2 Constructor & Destructor Documentation	399
10.103.3 Member Function Documentation	399
10.103.4 Event Documentation	412
10.104 CSerialPortNet Class Reference	412
10.104.1 Constructor & Destructor Documentation	412
10.104.2 Member Function Documentation	413
10.105 CStg200xBasicNet Class Reference	414
10.105.1 Detailed Description	418
10.105.2 Constructor & Destructor Documentation	418
10.105.3 Member Function Documentation	418
10.106 CStg200xDownloadBasicNet Class Reference	444
10.106.1 Detailed Description	446
10.106.2 Member Function Documentation	446
10.106.3 Property Documentation	453
10.107 CStg200xDownloadNet Class Reference	453
10.107.1 Detailed Description	454
10.107.2 Constructor & Destructor Documentation	454
10.107.3 Member Function Documentation	455
10.107.4 Event Documentation	459
10.108 CStimulusFunctionNet Class Reference	460
10.108.1 Constructor & Destructor Documentation	461
10.108.2 Member Function Documentation	462
10.108.3 Event Documentation	471
10.109 CSw2to64DeviceNet Class Reference	471
10.109.1 Detailed Description	472

10.109.2 Constructor & Destructor Documentation	472
10.109.3 Member Function Documentation	472
10.110 CTcxDeviceNet Class Reference	474
10.110.1 Detailed Description	476
10.110.2 Constructor & Destructor Documentation	476
10.110.3 Member Function Documentation	476
10.111 CTEERFunctionNet Class Reference	487
10.111.1 Detailed Description	489
10.111.2 Constructor & Destructor Documentation	489
10.111.3 Member Function Documentation	489
10.112 CTEERMachineDeviceNet Class Reference	497
10.112.1 Constructor & Destructor Documentation	497
10.112.2 Property Documentation	497
10.113 CUsbExceptionNet Class Reference	498
10.113.1 Detailed Description	498
10.113.2 Constructor & Destructor Documentation	498
10.113.3 Property Documentation	499
10.114 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference	499
10.114.1 Constructor & Destructor Documentation	499
10.114.2 Member Data Documentation	499
10.115 CW2100_FunctionNet Class Reference	499
10.115.1 Constructor & Destructor Documentation	501
10.115.2 Member Function Documentation	501
10.115.3 Property Documentation	506
10.116 CW2100_StimulatorFunctionNet Class Reference	507
10.116.1 Constructor & Destructor Documentation	508
10.116.2 Member Function Documentation	508
10.116.3 Member Data Documentation	513
10.116.4 Event Documentation	513
10.117 CW2100DacqGroupChannelSelectionNet Class Reference	513
10.117.1 Constructor & Destructor Documentation	514
10.118 CWarnerUssingDeviceNet Class Reference	514
10.118.1 Detailed Description	514
10.118.2 Constructor & Destructor Documentation	515
10.118.3 Property Documentation	515
10.119 CWarnerUssingFunctionNet Class Reference	515
10.119.1 Detailed Description	517
10.119.2 Constructor & Destructor Documentation	517
10.119.3 Member Function Documentation	518
10.120 CWarnerValveControllerDeviceNet Class Reference	532
10.120.1 Detailed Description	536
10.120.2 Constructor & Destructor Documentation	536

10.120.3 Member Function Documentation	536
10.120.4 Event Documentation	552
10.121 CWarnerValveControllerDeviceTesterFunctionNet Class Reference	555
10.121.1 Detailed Description	556
10.121.2 Constructor & Destructor Documentation	556
10.121.3 Member Function Documentation	556
10.122 CWClassicFunctionNet Class Reference	558
10.122.1 Constructor & Destructor Documentation	559
10.122.2 Member Function Documentation	559
10.123 CWirelessBaseFunctionNet Class Reference	563
10.123.1 Constructor & Destructor Documentation	563
10.123.2 Member Function Documentation	563
10.124 DeviceIdNet Struct Reference	563
10.124.1 Detailed Description	564
10.124.2 Constructor & Destructor Documentation	564
10.124.3 Member Function Documentation	564
10.124.4 Member Data Documentation	565
10.125 DigitalSource< digitalsourceenum > Class Template Reference	565
10.125.1 Constructor & Destructor Documentation	565
10.125.2 Member Function Documentation	565
10.125.3 Property Documentation	566
10.126 DriverVersionNet Class Reference	566
10.126.1 Detailed Description	567
10.126.2 Constructor & Destructor Documentation	567
10.126.3 Member Function Documentation	567
10.127 FirmwareDestinationNames Class Reference	571
10.127.1 Member Data Documentation	572
10.128 HeadStageIDType Class Reference	575
10.128.1 Member Enumeration Documentation	575
10.128.2 Constructor & Destructor Documentation	576
10.128.3 Member Function Documentation	576
10.128.4 Property Documentation	576
10.129 HeadstageIDTypeObject Class Reference	578
10.129.1 Constructor & Destructor Documentation	578
10.129.2 Member Function Documentation	578
10.129.3 Member Data Documentation	578
10.129.4 Property Documentation	579
10.130 HeadStageIDTypeState Class Reference	579
10.130.1 Property Documentation	579
10.131 mkfilterNet Class Reference	580
10.131.1 Member Function Documentation	580
10.132 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference	583

10.132.1 Member Function Documentation	584
10.133 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference	590
10.133.1 Member Function Documentation	590
10.134 CMeaAudioFunctionNet::s_setaudionet Struct Reference	590
10.134.1 Member Data Documentation	590
10.135 CStimulusFunctionNet::SidebandData Class Reference	590
10.135.1 Constructor & Destructor Documentation	591
10.135.2 Property Documentation	591
10.136 StgStatusNet Class Reference	591
10.136.1 Member Function Documentation	592
10.136.2 Member Data Documentation	592
10.137 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference	592
10.137.1 Constructor & Destructor Documentation	592
10.137.2 Property Documentation	593
10.138 usbSetupPacket_t Class Reference	593
10.138.1 Member Data Documentation	594
10.139 W2100_StimulusParametersNet Struct Reference	594
10.139.1 Member Data Documentation	594
Index	596

1 McsUsbNet.dll for MCS USB devices

1.1 Introduction

This DLL provides the .NET interface to MCS devices

- [STG200x & STG400x series stimulus generator](#)
- [MC_Card and MCS-USB-ME\(A\) data acquisition devices](#)
- [Other MCS USB devices](#)

1.2 System requirements

The DLL can be used with any .NET compatible language.

It requires the Microsoft C runtime libraries to be installed:

- MSVCR120.dll
- MSVCP120.dll

It also requires the USB driver and/or the MC_Card driver to be installed.

The simplest way to achieve this is to install the latest MC_Rack setup.

All examples assume that the [Mcs.Usb](#) namespace is loaded:

```
using namespace Mcs.Usb;
```

Include the file McsUsbNet.dll into the references of your project.

2 Other MCS USB devices

- For TCx device see [Mcs.Usb.CTcxDeviceNet](#)
- For SW2TO64 device see [MCS-USB-Sw2to64](#)
- For FluidControl device see [MCS FluidControl](#)

2.1 MCS-USB-Sw2to64 device

The class `UsbNetDll::CSw2to64DeviceNet` controls the setting of the switches in the MCS-USB-Sw2to64 device.

First construct an object of the class and connect with `UsbNetDll::CMcsUsbNet::Connect` to the device:

```
CSw2to64DeviceNet device = new CSw2to64DeviceNet();
device.Connect();
```

When the connect function is called without argument, the first MCS-USB-Sw2to64 found on the USB bus is used. When more than one MCS-USB-Sw2to64 is connected, this function can take for example the serial number of the required MCS-USB-Sw2to64 as argument.

The get the number of channels the device handles

```
int number = device.GetNumber();
```

Set all channel switches at once

```
byte z = 1;
byte[] pattern = new byte[number];
for(int i = 0; i < number; i++)
{
    pattern[i] = z; // pattern you want to switch this channel to
}
device.SetChannels(pattern);
```

Get all channel switches at once

```
byte[] pattern = device.GetChannels();
```

Set one channel switch

```
ushort index = 10;
byte pattern = 1;
device.SetChannel(index, pattern)
```

Get one channel switch

```
ushort index = 10;
byte pattern = device.GetChannel(index);
```

2.2 The MCS FluidControl Device

2.2.1 Introduction

The FluidControl Device can control up to 24 valves. The nominal voltage is 24V.

8 TTL level digital output ports are available and 8 TTL inputs can be read in.

The device has 8 ADC inputs with a range from 0V to 3.3V.

2.2.2 Access to the FluidControl device

The connection to a FluidControl device is established by `CFluidControlDevice::Connect`. When this function is called without argument, the first FluidControl device found on the USB bus is used. When more than one FluidControl devices are connected, this function can take the serial number of the required FluidControl device as argument.

```
CFluidControlDevice* m_dacq;
m_fluidcontrol = new CFluidControlDevice;
status = m_fluidcontrol->Connect();
```

The valves are controlled with the `CFluidControlDevice::SetValve` call. The argument given is a bit pattern of all valves which should be open.

The digital outputs can be controlled with the `CFluidControlDevice::SetDigout` call. Again, a bit pattern of all digital output pins which should be set to a logic high level is given as an argument.

The current state of the valves and the digital outputs can be read back with the `CFluidControlDevice::GetValve` and `CFluidControlDevice::GetDigout`

The command to read an ADC-Channel is `CFluidControlDevice::GetAdc`. Here the channelnummer which should be read in is given as an argument and the return value is the current Adc level.

The state of the digital inputs is read with the `CFluidControlDevice::GetDigin` call. Here the return value is the bit pattern of the digital inputs.

The connection to the device is closed with the `CFluidControlDevice::Disconnect` call.

3 MC_Card and MCS-USB-ME(A) Data Acquisition Devices (DAQ)

There are two different device types for MEA data aquisition (DAQ) devices. There are the MC_Card PCI device and the MCS-USB-ME(A) devices. Both types are supported by this library.

This library does not support the writing of the MCD file format!

The class `UsbNetDll::CMeaDeviceNet` is the base class for MEA devices and its implemantation as USB-MEA and MC_Card device classes.

In .NET both classes can be accessed by the contructor of the base class `UsbNetDll::CMeaDeviceNet` which constructs the correct underlying classes for the USB-MEA device (`UsbNetDll::CMeaUSBDeviceNet`) on the one hand or the MC_Card device (`UsbNetDll::CMC_CardDeviceNet`) on the other hand.

```
// MC_Card device with OnNewData callback
CMeaDeviceNet device = new CMeaDeviceNet (McsBusTypeEnumNet.MCS_PCI_BUS, null, OnNewData, OnError);
// or
// USB MEA device with OnChannelData callback
CMeaDeviceNet device = new CMeaDeviceNet (McsBusTypeEnumNet.MCS_PCI_BUS, OnChannelData, null, OnError);
```

For each type of device you can use any of the two types of data callbacks

A connection to a DAQ device is established by `UsbNetDll::CMcsUsbNet::Connect`. When this function is called without argument, the first DAQ device found on the USB bus is used. When more than one DAQ device of the specific type is connected, this function can take the serial number of the require DAQ device as an arument.

```
device.Connect();
```

Get the number of available analog hardware channels and set the number of channels to the maximum.

```
int hwchannels;
device.HWInfo().GetNumberOfHWADCCannels(out hwchannels);
device.SetNumberOfChannels(hwchannels);
int samplingrate = 1000; // MC_Card does not support all settings, please see MC_Rack for valid settings
device.SetSampleRate(samplingrate);
// Set the range according to the index (only valid for MC_Card)
device.SetVoltageRange(0);
```



```
device.EnableDigitalIn(true);
```

Get the layout to know how the data look like that you receive

```
int ana, digi, che, tim, block;
device.GetChannelLayout(out ana, out digi, out che, out tim, out block);
```

For the `UsbNetDll::OnChannelData` callback function you have to provide a definition of the channels you want to receive.

```
bool[] selChannels = new bool[block];
for (int i = 0; i < block; i++)
{
    selChannels[i] = true; // With true channel i is selected
    // selChannels[i] = false; // With false the channel i is deselected
}
channelblocksize = samplingrate / 10; // good choise for MC_Card
// queue size and threshold should be selected carefully
device.SetSelectedChannels(selChannels, 10 * channelblocksize, channelblocksize);
```

The `UsbNetDll::OnChannelData` callback function gets a callback for each channelblock that is defined. In this example a callback for each channel.

```
void OnChannelData(int CbHandle, int numSamples)
{
    int size_ret;
    ushort[] channeldata = device.ChannelBlock_ReadData(CbHandle, numSamples, out size_ret);
}
```

Or provide a `UsbNetDll::OnNewData` callback function. This callback gets the raw datastream as received

```
void OnNewData(CMcsUsbDacqDataPointer pointer, int numSamples)
{
    ushort[] data = device.ReadData(pointer, numSamples);
}
void OnError(String msg, int info)
{
    MessageBox.Show("Mea Device Error: " + msg);
}
```

see the `MeaExampleNet` in the example directory.

4 The MCS Robo Device

4.1 Introduction

Up to now two MCS devices exist that base on the Robo platform.

- The MCS Roboinject device is controlled by the `UsbNetDll::CRobolInjectDeviceNet` class.
- The MCS Roboocyte2 device is controlled by the `UsbNetDll::CRoboocyte2DeviceNet` class.

Both classes are derived from `UsbNetDll::CRoboDeviceNet`

5 STG200x & STG400x Series Stimulus Generators

5.1 Introduction

The STG200x & STG400x Series Stimulus Generators have two distinct modes of operation, the [Download mode](#) and the [Streaming mode](#).

5.2 Download mode

The Download mode is the "classic" mode of operation, as used by the MC Stimulus software. In this mode, one or multiple waveforms are defined in PC memory and downloaded to the STG. The waveforms are stored in STG device onboard memory and can be sent to the analog and sync outputs once or multiple times. The STG can operate independantly from the PC (without computer connection) after the download. Output is triggered either by the front panel start/stop button, the digital trigger inputs or under software control.

In the Download mode, there are up to eight independent triggers available (depending on the device). The user can assign each of the analog outputs and sync (digital) outputs to any of the triggers.

The analog output waveform is stored sample by sample in the STG memory. To reduce memory usage, this data can be compressed: whenever a given output value is to be held for more than one sample period, it has only to be given once. The user can define the number of sample periods for that a pattern should remain active. Compression is done for each channel independantly of the others, thus the algorithm to compress the data is very easy to implement.

A new feature of the Download mode is the segmentation of the STG memory. The onboard memory can be devided into up to 100 segments. Each segment can hold its own waveform pattern. Under software control, the user can switch between the defined segments within milliseconds. Another option is to use the four trigger inputs to select between four predefined segments. This option is accessible from the MC_Stimulus Software as the "Multi-File mode", and can start each of up to four defined waveforms within microseconds. This feature allows a predefined flexible response (feedback) to recorded data.

UsbNetDll::CStg200xDownloadNet is the class for using the STG in download mode.

5.2.1 Memory Layout and Trigger Setup

The class to be used for the Download mode is UsbNetDll::CStg200xDownloadNet, which is derived from UsbNetDll::CStg200xBasicNet. You can add a poll handler delegate (UsbNetDll::OnStg200xPollStatus) to the constructor UsbNetDll::CStg200xDownloadNet::CStg200xDownloadNet.

The connection to a STG is established by UsbNetDll::CMcsUsbNet::Connect. When this function is called without argument, the first STG found on the USB bus is used. When more than one STG is connected, this function can take for example the serial number of the required STG as argument.

```
CStg200xDownloadNet device = new CStg200xDownloadNet();
device.Connect();
```

To use the Download mode, the memory layout of the STG200x can be set up, if the default is not sufficient. The total amount of memory available in the STG is obtained by the UsbNetDll::CStg200xDownloadNet::GetTotalMemory call. With UsbNetDll::CStg200xDownloadNet::SendSegmentDefine the segment sizes are assigned.

```
uint32_t memory = device.GetTotalMemory(); // obtain total memory available
uint[] segmentmemory = new uint[2];      // each segments has half of total memory
segmentmemory[0] = memory / 2;
segmentmemory[1] = memory / 2;
device.SendSegmentDefine(segmentmemory); // setup the STG
```

Next, for each segment, one has to assign the amount of memory to be used for each channel and sync output. This is done by UsbNetDll::CStg200xDownloadBasicNet::SetCapacity. Its arguments contain a list of memory sizes, with one entry per channel and one entry per sync output. Again, the total memory assigned to the channels and sync outputs must not exceed the memory assigned to the segment.

```
uint32_t nchannels = device.GetNumberOfAnalogChannels();
uint32_t nsync = device.GetNumberOfSyncoutChannels();
uint[] channel_cap = new uint[nchannels];
uint[] syncout_cap = new uint[nsync];
for (int i = 0; i < 2; i++) // for each segment
{
    device.SendSegmentSelect((uint32_t)i); // switch to segment
    uint32_t segment_mem = device.GetMemory(); // get memory available in this segment
    for(int j = 0; j < nchannels; j++)
    {
        channel_cap[j] = segment_mem/(nchannels+nsync); // devide memory amount to all channels
    }
}
```

```

for(int j = 0; j < nsync; j++)
{
    syncout_cap[j] = segment_mem/(nchannels+nsync); // and all sync outs.
}
device.SetCapacity(channel_cap, syncout_cap); // define memory for current segment
}

```

Before the STG can start, the trigger has to be configured. This is done by the `UsbNetDll::CStg200xDownloadNet::SetupTrigger` call. Its arguments are a list of channelmaps, syncoutmaps and repeats, one for each of the four available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and syncout 1 to trigger 1 and channel 3 to trigger 2 use:

```

uint32_t TriggerInputs = device.GetNumberOfTriggerInputs();
uint[] channelmap = new uint[TriggerInputs];
uint[] syncoutmap = new uint[TriggerInputs];
uint[] repeat = new uint[TriggerInputs];
for (int i = 0; i < TriggerInputs; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    repeat[i] = 0;
}
// Trigger 0
channelmap[0] = 1; // Channel 1
syncoutmap[0] = 1; // Syncout 1
repeat[0] = 0; // forever
// Trigger 1
channelmap[1] = 4; // Channel 3
device.SetupTrigger(channelmap, syncoutmap, repeat);

```

For the STG400x series you have to set the output mode of the channels. `UsbNetDll::CStg200xDownloadNet::SetVoltageMode` interprets the values as voltages. `UsbNetDll::CStg200xDownloadNet::SetCurrentMode` as currents.

```

// Only meaningful for STG400x
device.SetVoltageMode();

```

For each segment, data can be sent to each of the defined channels and sync outputs using the `UsbNetDll::CStg200xDownloadNet::SendChannelData` and `UsbNetDll::CStg200xDownloadNet::SendSyncData` calls. channeldata and syncdata are a list of analog and digital samples as a list of two byte values (unsigned short). Multiple calls to `UsbNetDll::CStg200xDownloadNet::SendChannelData` and `UsbNetDll::CStg200xDownloadNet::SendSyncData` to the same channel append data to that channel.

If the Multi-File mode of the STG is enabled using the `UsbNetDll::CStg200xDownloadNet::EnableMultiFileMode` call, the four trigger inputs are used to switch between four segments. A hardware trigger signal (TTL) on trigger input 1 selects the first segment and starts all pulses in this segment. Thus with the Multi-File mode, one can predefined four stimulus patterns and switch between them without a connection to the PC.

The STG200x series has an analog resolution of 13 bits, thus the analog data contains the information in bits 0 to 12 of each sample. Bits 13 to 15 have to be 0.

```

int DACResolution = device.GetDACResolution();
// Data for Channel 0
{
    device.ClearChannelData(0);
    double factor = 0.1;
    const int l = 1000;
    ushort[] pData = new ushort[l];
    UInt64_t[] tData = new UInt64_t[l];
    for (int i = 0; i < l; i++)
    {
        // calculate Sin-Wave
        double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
            Math.Sin(2.0 * (double)i * Math.PI / (double)l);
        // calculate sign
        pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
            (int)Math.Pow(2, DACResolution - 1));
        tData[i] = (UInt64_t)20; // duration in µs
    }
    device.SendChannelData(0, pData, tData);
}
// Data for Channel 3
{
    device.ClearChannelData(2);
    double factor = 0.1;
    const int l = 700;
    // without compression
    ushort[] pData = new ushort[l];
    UInt64_t[] tData = new UInt64_t[l];
    for (int i = 0; i < l; i++)

```

```

{
    // calculate Sin-Wave
    double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
        Math.Sin(2.0 * (double)i * Math.PI / (double)1);
    // calculate sign
    pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
        (int)Math.Pow(2, DACResolution - 1));
    tData[i] = (UInt64_t)20; // duration in µs
}
device.SendChannelData(2, pData, tData);
}
// Data for Sync 0
{
    device.ClearSyncData(0);
    ushort[] pData = new ushort[1000];
    UInt64_t[] tData = new UInt64_t[1000];
    for (int i = 0; i < 1000; i++)
    {
        pData[i] = (ushort)(i&1);
        tData[i] = 20;
    }
    device.SendSyncData(0, pData, tData);
}

```

Start the trigger by pushing the front button or by software

```

// Start Trigger 1 and 2
device.SendStart(1 + 2); // Trigger 1 und 2

```

see the StgDownloadExampleNet in the example directory.

5.3 Streaming mode

The other mode of operation is the Streaming mode. Here the analog output is sent to the STG device in "real time". The PC has to be connected to the STG all the time. The data that is sent to the analog output is downloaded from the PC to the STG on the fly.

The Streaming mode is useful for applications where flexible feedback is needed as well for applications where very long waveforms which are not repeated (such as white noise) are used.

The Streaming mode works by use of two ring buffers which hold data. One is in PC memory and managed by the DLL, and one is in on-board STG memory. Data is transferred from PC memory to the STG via the USB bus in time slices of one millisecond.

The user can define both the size of the ring buffer in DLL memory and in the STG memory. Once the Streaming mode is started, the STG request data from the PC. The data rate from PC to STG is variable and controlled by the STG. The STG request data from the PC at a rate to keep its internal ringbuffer at about half full.

It is the responsibility of the user to keep the ring buffer in the memory of the PC filled, so the DLL can supply sufficient data to the STG. To do so, the Windows DLL allows to define a "callback" function which is called whenever new data is needed, or more precise, as soon as the ring buffer in the memory of the PC falls below the user defined threshold.

Small buffers have the advantage of a low latency between data generation in the callback function and its output as a analog signal from the STG. However for low latency to work, the user-written callback function has to be fast and to produce a steady flow of data.

In the Streaming mode, all triggers are available as well. Each of the eight analog and sync outputs can be assigned to one of the triggers.

The output rate is user defined with a maximum of 50 kHz

UsbNetDll::CStg200xStreamingNet is the class for using the STG in streaming mode.

5.3.1 Memory Layout and Trigger Setup

With the constructor for `UsbNetDll::CStg200xStreamingNet::CStg200xStreamingNet`, the name of the callback function for the data handler is provided. The data handler function is called automatically, whenever the STG needs new data. This data is first written to a ring buffer in the memory of the PC. The size for this ring buffer is defined as first argument in the constructor. The user provided delegate gets the trigger number which needs new data as argument

```
CStg200xStreamingNet device = new CStg200xStreamingNet(10000, dataHandler, errorHandler);
```

The callback function, which is defined in the constructor, is called whenever the STG needs new data for a trigger, or more precise, whenever the ring buffer in PC memory falls below the defined threshold.

The user can query the amount of space available for queuing by use of the `UsbNetDll::CStg200xStreamingNet::GetDataQueueSpace` call. Its return value is the number of samples that can be send to the STG.

User code is required to fill an array analog and sync out data, sample by sample for up to the maximum number of samples as obtained by `UsbNetDll::CStg200xStreamingNet::GetDataQueueSpace` or `UsbNetDll::CStg200xStreamingNet::GetSyncoutQueueSpace`.

The values for the analog outputs are 16 bits signed integers. The lower bits are truncated according to the resolution of the STG. This behaviour is different to the behaviour in [download mode](#).

Note: Compression as described in the [download mode](#) can NOT be used for the streaming mode.

The new data is sent to the STG by using the `UsbNetDll::CStg200xStreamingNet::EnqueueData` call.

```
void dataHandler(uint32_t trigger)
{
    double factor = 1;
    if (trigger == 0) // Callback for Trigger 1
    {
        // Handle Channel 1
        uint32_t channel = 0;
        for ( ; ; )
        {
            uint32_t space = device.GetDataQueueSpace(channel);
            if (space < 1000)
                break;
            short[] data = new short[1000];
            for (int i = 0; i < 1000; i++)
            {
                // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
                double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) *
                    Math.Sin(2.0 * (double)i * Math.PI / (double)1000);
                data[i] = (short)sin;
            }
            uint32_t enqueued = device.EnqueueData(channel, data);
        }
    }
    // Handle Channel 3
    uint32_t channel = 2;
    for ( ; ; )
    {
        uint32_t space = device.GetDataQueueSpace(channel);
        if (space < 700)
            break;
        short[] data = new short[700];
        for (int i = 0; i < 700; i++)
        {
            // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
            double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) *
                Math.Sin(2.0 * (double)i * Math.PI / (double)700);
            data[i] = (short)sin;
        }
        uint32_t enqueued = device.EnqueueData(channel, data);
    }
}
// Handle Syncout 1
uint32_t channel = 0;
for ( ; ; )
{
    uint32_t space = device.GetSyncoutQueueSpace(channel);
    if (space < 1000)
        break;
    ushort[] data = new ushort[1000];
    for (int i = 0; i < 1000; i++)
```

```

        data[i] = (ushort)(i & 1);
        uint32_t enqueued = device.EnqueueSyncout(channel, data);
    }
}
}
}
void errorHandler()
{
}

```

The connection to a STG is established by `UsbNetDll::CMcsUsbNet::Connect`. When this function is called without argument, the first STG found on the USB bus is used. When more than one STG is connected, this function can take the serial number of the required STG as an argument.

```
device.Connect();
```

With enabling or disabling the continuous mode it can be selected how the STG handles an "out of data" situation.

When `UsbNetDll::CStg200xStreamingNet::EnableContinuousMode` is used, the STG does not stop when it runs out of data, but it keeps running and sends a zero voltage to its outputs.

When `UsbNetDll::CStg200xStreamingNet::DisableContinuousMode` is used, the STG stops when it runs out of data. It has to be retriggered to resume the output.

```
device.EnableContinuousMode();
```

`UsbNetDll::CStg200xStreamingNet::SetOutputRate` is used to set the sampling rate.

```
device.SetOutputRate(50000);
```

To use the Streaming mode, the memory layout of the STG has to be set up. To total amount of memory available in the STG is obtained by the `UsbNetDll::CStg200xStreamingNet::GetTotalMemory` call.

This memory can be assigned to four ring buffers (one per trigger) which buffer the data received from the PC via USB cable. This is done with the `CStg200xStreaming::SetCapacity` call. The total amount of memory must not exceed the total memory size as obtained by `UsbNetDll::CStg200xStreamingNet::GetTotalMemory`.

This internal ring buffer is crucial for proper operation of the Streaming mode. The size of the ring buffer determines the latency of the Streaming mode. The firmware of the STG requests data from the PC in order to keep the ring buffer about half full. Thus the average latency is:

$$\text{latency} = (\text{ringbuffersize in bytes} / 4) / \text{output rate}$$

If the ring buffer size is too big, the latency of the STG might be too long. If the ring buffer size is too low, an overflow or underflow of data in the STG ringbuffer might occur, resulting in data jumps of the output signals or the "out of data" situation described earlier.

The following example divides the total memory equally among the four triggers:

```

uint32_t dwMemory = device.GetTotalMemory(); // obtain total memory available
uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] stg_triggercapacity = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
    stg_triggercapacity[i] = dwMemory / ntrigger;
device.SetCapacity(stg_triggercapacity); // setup the STG

```

or fixed memory sizes:

```

uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] stg_triggercapacity = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
    stg_triggercapacity[i] = 50000;
device.SetCapacity(stg_triggercapacity);

```

Before the STG can start, the trigger has to be configured. This is done by the `UsbNetDll::CStg200xStreamingNet::SetupTrigger` call. Its arguments are a list of channelmaps, syncoutmaps, digoutmap, autostart and callback, with one entry for each of the available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and 3 and syncout 1 to trigger 1 use:

```

uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] channelmap = new uint[ntrigger];
uint[] syncoutmap = new uint[ntrigger];
uint[] digoutmap = new uint[ntrigger];
uint[] autostart = new uint[ntrigger];

```

```
uint[] callback_threshold = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    digoutmap[i] = 0;
    autostart[i] = 0;
    callback_threshold[i] = 0;
}
channelmap[0] = 0x1 + 0x4; // Channel 1 und Channel 3 to Trigger 1
syncoutmap[0] = 0x1; // Syncout 1 to Trigger 1
autostart[0] = 1;
callback_threshold[0] = 50; // 50% of buffer size
device.SetupTrigger(channelmap, syncoutmap, digoutmap, autostart, callback_threshold);
device.StartLoop();
System.Threading.Thread.Sleep(1000); // Give StartLoop some time
```

Start Trigger by pushing the front button or by Software
device.SendStart(1);

see the StgStreamingExampleNet in the example directory.

6 Namespace Index

6.1 Namespace List

Here is a list of all namespaces with brief descriptions:

Mcs	21
Mcs::Usb	21

7 Hierarchical Index

7.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

CW2100_FunctionNet::AudioChannelsNet	27
BatteryState	28
CCreateFilterNet	46
BesselFilterHighPassNet	28
BesselFilterLowPassNet	29
ButterworthFilterHighPassNet	29
ButterworthFilterLowPassNet	30
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet >	55
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumNet >	55
CDeviceGroupChannelInfoNet	54
CDeviceGroupChannelInfoTemplateNet< int >	55

CDeviceGroupChannelInfoGenericNet	53
CDeviceGroupChannelInfoTemplateNet< SCUDacqGroupChannelEnumNet >	55
CDeviceGroupChannelInfoSCUNet	54
CDeviceGroupChannelInfoTemplateNet< W2100DacqGroupChannelEnumNet >	55
CDeviceGroupChannelInfoW2100Net	56
CFilterCoefficientsNet	63
CFilterPropertyNet	70
CMcsUsbDacqNet::CHWInfo	101
CMcsUsbFunctionNet	215
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumNet, DacqGroupChannelEnum, CDeviceGroupChannelInfoNet >	51
CDacqGroupChannelSelectionNet	50
CDacqGroupChannelSelectionTemplateNet< int, int, CDeviceGroupChannelInfoGenericNet >	51
CDacqGroupChannelGenericSelectionNet	49
CDacqGroupChannelSelectionTemplateNet< SCUDacqGroupChannelEnumNet, SCUDacqGroupChannelEnum, CDeviceGroupChannelInfoSCUNet >	51
CSCUDacqGroupChannelSelectionNet	396
CDacqGroupChannelSelectionTemplateNet< W2100DacqGroupChannelEnumNet, W2100DacqGroupChannelEnum, CDeviceGroupChannelInfoW2100Net >	51
CW2100DacqGroupChannelSelectionNet	513
CCMOSMea_FunctionNet	32
CDacCalibrationFunctionNet	48
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >	51
CDigOutStimulatorFunctionNet	56
CFilterConfigurationNet	67
CFilterConfigurationRegisterNet	68
CIntanMea_FunctionNet	103
CInterfaceboardFunctionNet	105
CMcsBus_AxisParametersNet	114
CMcsBus_ExtensionNet	116
CMcsBus_FYIExtensionNet	117
CMcsBus_MotorControlNet	119
CMcsBus_SensorNet	136

CMcsBus_TempSensorNet	146
CMcsBus_VoltageModeNet	148
CMcsBusNet	153
CMcsUsbDeviceStatePushFunctionNet	205
CMultiwellCallbackFunctionNet	288
CSCUFunctionNet	397
CMEA2100x256FunctionNet	245
CMeaAudioFunctionNet	247
CMeaDigitalDataFunctionNet	268
CMeaFeedbackFunctionNet	270
CMeFunctionNet	279
CMultiwellOptoStimFunctionNet	297
CPPCFunctionNet	324
CPPS_FunctionNet	331
CProgramPressureCurveNet	336
CPulseGeneratorFunctionNet	338
CRFFunctionNet	346
CRobo_FYIProgram_FunctionNet	350
CRobo_FYITemp_FunctionNet	352
CStimulusFunctionNet	460
CTEERFunctionNet	487
CW2100_StimulatorFunctionNet	507
CWarnerUssingFunctionNet	515
CWarnerValveControllerDeviceTesterFunctionNet	555
CWirelessBaseFunctionNet	563
CW2100_FunctionNet	499
CWClassicFunctionNet	558
CMcsUsbFunctionPointerContainer	217
CMcsUsbListEntryNet	217
CMcsUsbListNet	221
CMcsUsbNet	224
CExternDTesterDeviceNet	61

CFluidControlDeviceNet	72
CGenericDevelopDeviceNet	79
CGilsonDeviceNet	96
CMcsUsbDacqNet	157
CMeaDeviceNet	260
CMeaUSBDeviceNet	278
CCMOSMeaDeviceNet	43
CHLADacqNet	99
CLIH3DeviceNet	107
CMultiwellDeviceNet	291
CWarnerUssingDeviceNet	514
COctoPotDeviceNet	302
CRoboDacqNet	354
CMcsUsbDeviceStatePushNet	206
CWarnerValveControllerDeviceNet	532
CMcsUsbFactoryNet	207
CMeaCleanDeviceNet	250
CMeaCoatDeviceNet	254
CMeaImpedanceDeviceNet	274
CMeaSwitchDeviceNet	276
CChannelTestDeviceNet	31
CMultiBatteryChargerDeviceNet	281
CNF_GenDeviceNet	302
COkuvisionStimulatorDeviceNet	306
CPathIdentDeviceNet	311
CPedoterDeviceNet	312
CPeristalticPumpDeviceNet	313
CPgaDeviceNet	314
CPositionIIDeviceNet	316
CPositionImpDeviceNet	319
CPPCDeviceNet	322
CPPS_DeviceNet	330

CRadioControlledDevicesNet	341
CRetinaLedDeviceNet	344
CRoboDeviceNet	367
CEncapsulatorDeviceNet	60
CFYIDeviceNet	78
CHiClampDeviceNet	98
CHLADeviceNet	100
CMeasureTableDeviceNet	275
CPatchServerDeviceNet	310
CPPSDeviceNet	336
CRoboStatorDeviceNet	387
CRoboInjectDeviceNet	385
CRoboocyte2DeviceNet	385
CTEERMachineDeviceNet	497
CRoboFluidDeviceNet	381
CSafeISDeviceNet	393
CSerialPortNet	412
CStg200xBasicNet	414
CStg200xDownloadBasicNet	444
CStg200xDownloadNet	453
CSw2to64DeviceNet	471
CTcxDeviceNet	474
CMcsUsbPointerContainer	245
CCMOSMeaDeviceNet::CRegionOfInterestRect	343
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	499
DeviceIdNet	563
DigitalSource< digitalsourceenum >	565
DriverVersionNet	566
Exception	
CUsbExceptionNet	498
FirmwareDestinationNames	571
HeadstageIDTypeObject	578

HeadStageIDTypeState	579
Comparable	
HeadStageIDType	575
mkfilterNet	580
CRoboDeviceNet::RoboMainLowLevelCommands	583
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	590
CMeaAudioFunctionNet::s_setaudionet	590
CStimulusFunctionNet::SidebandData	590
StgStatusNet	591
stgstreaming	
CStg200xBasicNet	414
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	592
usbSetupPacket_t	593
W2100_StimulusParametersNet	594

8 Class Index

8.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

CW2100_FunctionNet::AudioChannelsNet	27
BatteryState	28
BesselFilterHighPassNet	28
BesselFilterLowPassNet	29
ButterworthFilterHighPassNet	29
ButterworthFilterLowPassNet	30
CChannelTestDeviceNet	31
CCMOSMea_FunctionNet	32
CCMOSMeaDeviceNet	43
CCreateFilterNet	46
CDacCalibrationFunctionNet	
48	
CDacqGroupChannelGenericSelectionNet	49
CDacqGroupChannelSelectionNet	50

CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplateNet, 51	
CDeviceGroupChannelInfoGenericNet	53
CDeviceGroupChannelInfoNet	54
CDeviceGroupChannelInfoSCUNet	54
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet >	55
CDeviceGroupChannelInfoW2100Net	56
CDigOutStimulatorFunctionNet	
CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class	56
CEncapsulatorDeviceNet	
CEncapsulatorDeviceNet is the to control the MCS HiClamp device	60
CExternDTesterDeviceNet	
CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)	61
CFilterCoefficientsNet	63
CFilterConfigurationNet	67
CFilterConfigurationRegisterNet	68
CFilterPropertyNet	70
CFluidControlDeviceNet	
CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device	72
CFYIDeviceNet	
CFYIDeviceNet is the class to control the MCS FYI device	78
CGenericDevelopDeviceNet	
CGenericDevelopDeviceNet is the class to use during development of a new device	79
CGilsonDeviceNet	
CGilsonDeviceNet is the class to control a Gilson device	96
CHiClampDeviceNet	
CHiClampDeviceNet is the to control the MCS HiClamp device	98
CHLADacqNet	99
CHLADeviceNet	
CHLADeviceNet is the to control the MCS HLA device	100
CMcsUsbDacqNet::CHWInfo	
Class to provide hardware information about the device	101
CIntanMea_FunctionNet	103
CInterfaceboardFunctionNet	
CInterfaceboardFunctionNet is the class to control the Interfaceboard	105
CLIH3DeviceNet	
CLIH3DeviceNet is the class to access the HEKA LIH3 device	107
CMcsBus_AxisParametersNet	114

CMcsBus_ExtensionNet	116
CMcsBus_FYIExtensionNet	117
CMcsBus_MotorControlNet	119
CMcsBus_SensorNet	136
CMcsBus_TempSensorNet	146
CMcsBus_VoltageModeNet	148
CMcsBusNet	153
CMcsUsbDacqNet	
Base class for data acquisition devices	157
CMcsUsbDeviceStatePushFunctionNet	205
CMcsUsbDeviceStatePushNet	206
CMcsUsbFactoryNet	207
CMcsUsbFunctionNet	215
CMcsUsbFunctionPointerContainer	217
CMcsUsbListEntryNet	
McsUsbListEntryNet identifies a connected device	217
CMcsUsbListNet	
Class to handle a list of connected MCS USB devices	221
CMcsUsbNet	
Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class	224
CMcsUsbPointerContainer	245
CMEA2100x256FunctionNet	
CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "↵ Stg200xNet.h" to resolve documentation reference	245
CMeaAudioFunctionNet	247
CMeaCleanDeviceNet	
CMeaCleanDeviceNet is the class to access the MEA Clean device	250
CMeaCoatDeviceNet	
CMeaCoatDeviceNet is the class to access the MEA Coat device	254
CMeaDeviceNet	
Base class for MEA data acquisition devices	260
CMeaDigitalDataFunctionNet	268
CMeaFeedbackFunctionNet	270
CMeaImpedanceDeviceNet	274
CMeasureTableDeviceNet	
CMeasureTableDeviceNet is the to control the MCS HLA device	275

CMeaSwitchDeviceNet	
The class to control the USB-MEA-Switch	276
CMeaUSBDeviceNet	
Class for data acquisition via ME and MEA USB amplifiers	278
CMeFunctionNet	
279	
CMultiBatteryChargerDeviceNet	
CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device	281
CMultiwellCallbackFunctionNet	
CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator	288
CMultiwellDeviceNet	
CMultiwellDeviceNet is the class to access the HEKA LIH3 device	291
CMultiwellOptoStimFunctionNet	
CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device	297
CNF_GenDeviceNet	302
COctoPotDeviceNet	302
CokuvisionStimulatorDeviceNet	306
CPatchServerDeviceNet	
CPatchServerDeviceNet is the class to control the MCS PatchServer device	310
CPathIdentDeviceNet	311
CPedoterDeviceNet	
312	
CPeristalticPumpDeviceNet	
CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump	313
CPgaDeviceNet	314
CPositionIIDeviceNet	
CPositionIIDeviceNet is the class to control PositionII devices	316
CPositionImpDeviceNet	
CPositionImpDeviceNet is the class to access the Position/Imp devices	319
CPPCDeviceNet	322
CPPCFunctionNet	
CPPCFunctionNet is the class to access the PPC (high precision Patch Peristaltic patch Pump	324
CPPS_DeviceNet	330
CPPS_FunctionNet	331
CPPSDeviceNet	
CPPS4plus1DeviceNet is the to control the MCS HLA device	336
CProgramPressureCurveNet	
CProgramPressureCurveNet is the class to program pressure curves	336

CPulseGeneratorFunctionNet	
CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking	338
CRadioControlledDevicesNet	341
CCMOSMeaDeviceNet::CRegionOfInterestRect	343
CRetinaLedDeviceNet	344
CRFFFunctionNet	
CRFFFunctionNet is the class to control RF devices	346
CRobo_FYIProgram_FunctionNet	350
CRobo_FYITemp_FunctionNet	352
CRoboDacqNet	354
CRoboDeviceNet	
CRoboDeviceNet is the base class for all Robo platform based devices	367
CRoboFluidDeviceNet	381
CRoboInjectDeviceNet	
CRoboInjectDeviceNet is the to control the MCS RoboInject device	385
CRoboocyte2DeviceNet	
CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device	385
CRoboStatorDeviceNet	387
CSafeISDeviceNet	
393	
CSCUDacqGroupChannelSelectionNet	396
CSCUFunctionNet	
CSCUFunctionNet is the class to control the SCU device	397
CSerialPortNet	412
CStg200xBasicNet	
Base class for the Stg200x	414
CStg200xDownloadBasicNet	
CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device	444
CStg200xDownloadNet	
Main class for the STG download mode This class implements the STG download mode interface.	453
CStimulusFunctionNet	460
CSw2to64DeviceNet	
The class to control the MCS-USB-Sw2to64 device	471
CTcxDeviceNet	
Class to control a Temperature Controller (TCX)	474
CTEERFunctionNet	
CTEERFunctionNet is the class to control the TEER device	487

CTEERMachineDeviceNet	497
CUsbExceptionNet Exception class that is thrown in case of an USB error	498
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	499
CW2100_FunctionNet	499
CW2100_StimulatorFunctionNet	507
CW2100DacqGroupChannelSelectionNet	513
CWarnerUssingDeviceNet CWarnerUssingDeviceNet is the class to control the Ussing device	514
CWarnerUssingFunctionNet CWarnerUssingFunctionNet is the class to control the Ussing device	515
CWarnerValveControllerDeviceNet CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller	532
CWarnerValveControllerDeviceTesterFunctionNet CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester	555
CWClassicFunctionNet	558
CWirelessBaseFunctionNet	563
DeviceIdNet Device Id	563
DigitalSource< digitalsourceenum >	565
DriverVersionNet Class gives firmware versions of the device's firmware destinations	566
FirmwareDestinationNames	571
HeadStageIDType	575
HeadstageIDTypeObject	578
HeadStageIDTypeState	579
mkfilterNet	580
CRoboDeviceNet::RoboMainLowLevelCommands	583
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	590
CMeaAudioFunctionNet::s_setaudionet	590
CStimulusFunctionNet::SidebandData	590
StgStatusNet	591
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	592
usbSetupPacket_t	593
W2100_StimulusParametersNet	594

9 Namespace Documentation

9.1 Mcs Namespace Reference

Namespaces

- [Usb](#)

9.2 Mcs::Usb Namespace Reference

Classes

- class [BatteryState](#)
- class [BesselFilterHighPassNet](#)
- class [BesselFilterLowPassNet](#)
- class [ButterworthFilterHighPassNet](#)
- class [ButterworthFilterLowPassNet](#)
- class [CChannelTestDeviceNet](#)
- class [CCMOSMea_FunctionNet](#)
- class [CCMOSMeaDeviceNet](#)
- class [CCreateFilterNet](#)
- class [CDacCalibrationFunctionNet](#)
- class [CDacqGroupChannelGenericSelectionNet](#)
- class [CDacqGroupChannelSelectionNet](#)
- class [CDacqGroupChannelSelectionTemplateNet](#)
- class [CDeviceGroupChannelInfoGenericNet](#)
- class [CDeviceGroupChannelInfoNet](#)
- class [CDeviceGroupChannelInfoSCUNet](#)
- class [CDeviceGroupChannelInfoTemplateNet](#)
- class [CDeviceGroupChannelInfoW2100Net](#)
- class [CDigOutStimulatorFunctionNet](#)
[CDigOutStimulatorFunctionNet](#) is the class of the DigOut stimulator function class.
- class [CEncapsulatorDeviceNet](#)
[CEncapsulatorDeviceNet](#) is the to control the MCS HiClamp device
- class [CExternDTesterDeviceNet](#)
[CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)
- class [CFilterCoefficientsNet](#)
- class [CFilterConfigurationNet](#)
- class [CFilterConfigurationRegisterNet](#)
- class [CFilterPropertyNet](#)
- class [CFluidControlDeviceNet](#)
[CFluidControlDeviceNet](#) is the class to control MCS FluidControl (FCB and FCX) device.
- class [CFYIDeviceNet](#)
[CFYIDeviceNet](#) is the class to control the MCS FYI device
- class [CGenericDevelopDeviceNet](#)
[CGenericDevelopDeviceNet](#) is the class to use during development of a new device.
- class [CGilsonDeviceNet](#)
[CGilsonDeviceNet](#) is the class to control a Gilson device.
- class [CHiClampDeviceNet](#)
[CHiClampDeviceNet](#) is the to control the MCS HiClamp device
- class [CHLADacqNet](#)

- class [CHLADeviceNet](#)
CHLADeviceNet is the to control the MCS HLA device
- class [CIntanMea_FunctionNet](#)
- class [CInterfaceboardFunctionNet](#)
CInterfaceboardFunctionNet is the class to control the Interfaceboard
- class [CLIH3DeviceNet](#)
CLIH3DeviceNet is the class to access the HEKA LIH3 device.
- class [CMcsBus_AxisParametersNet](#)
- class [CMcsBus_ExtensionNet](#)
- class [CMcsBus_FYIExtensionNet](#)
- class [CMcsBus_MotorControlNet](#)
- class [CMcsBus_SensorNet](#)
- class [CMcsBus_TempSensorNet](#)
- class [CMcsBus_VoltageModeNet](#)
- class [CMcsBusNet](#)
- class [CMcsUsbDacqNet](#)
Base class for data acquisition devices.
- class [CMcsUsbDeviceStatePushFunctionNet](#)
- class [CMcsUsbDeviceStatePushNet](#)
- class [CMcsUsbFactoryNet](#)
- class [CMcsUsbFunctionNet](#)
- class [CMcsUsbFunctionPointerContainer](#)
- class [CMcsUsbListEntryNet](#)
McsUsbListEntryNet identifies a connected device.
- class [CMcsUsbListNet](#)
Class to handle a list of connected MCS USB devices.
- class [CMcsUsbNet](#)
Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.
- class [CMcsUsbPointerContainer](#)
- class [CMEA2100x256FunctionNet](#)
CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference
- class [CMeaAudioFunctionNet](#)
- class [CMeaCleanDeviceNet](#)
CMeaCleanDeviceNet is the class to access the MEA Clean device.
- class [CMeaCoatDeviceNet](#)
CMeaCoatDeviceNet is the class to access the MEA Coat device.
- class [CMeaDeviceNet](#)
Base class for MEA data acquisition devices.
- class [CMeaDigitalDataFunctionNet](#)
- class [CMeaFeedbackFunctionNet](#)
- class [CMeaImpedanceDeviceNet](#)
- class [CMeasureTableDeviceNet](#)
CMeasureTableDeviceNet is the to control the MCS HLA device
- class [CMeaSwitchDeviceNet](#)
The class to control the USB-MEA-Switch.
- class [CMeaUSBDeviceNet](#)
Class for data acquisition via ME and MEA USB amplifiers
- class [CMeFunctionNet](#)
- class [CMultiBatteryChargerDeviceNet](#)
CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

- class [CMultiwellCallbackFunctionNet](#)
CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator
- class [CMultiwellDeviceNet](#)
CMultiwellDeviceNet is the class to access the HEKA LIH3 device.
- class [CMultiwellOptoStimFunctionNet](#)
CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device
- class [CNF_GenDeviceNet](#)
- class [COctoPotDeviceNet](#)
- class [COKuvisionStimulatorDeviceNet](#)
- class [CPatchServerDeviceNet](#)
CPatchServerDeviceNet is the class to control the MCS PatchServer device
- class [CPathIdentDeviceNet](#)
- class [CPedoterDeviceNet](#)
- class [CPeristalticPumpDeviceNet](#)
CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.
- class [CPgaDeviceNet](#)
- class [CPositionIIDeviceNet](#)
CPositionIIDeviceNet is the class to control PositionII devices
- class [CPositionImpDeviceNet](#)
CPositionImpDeviceNet is the class to access the Position/Imp devices
- class [CPPCDeviceNet](#)
- class [CPPCFunctionNet](#)
CPPCFunctionNet is the class to access the PPC (high precision Patch Peristaltic patch Pump
- class [CPPS_DeviceNet](#)
- class [CPPS_FunctionNet](#)
- class [CPPSDeviceNet](#)
CPPS4plus1DeviceNet is the to control the MCS HLA device
- class [CProgramPressureCurveNet](#)
CProgramPressureCurveNet is the class to program pressure curves
- class [CPulseGeneratorFunctionNet](#)
CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking
- class [CRadioControlledDevicesNet](#)
- class [CRetinaLedDeviceNet](#)
- class [CRFFFunctionNet](#)
CRFFFunctionNet is the class to control RF devices
- class [CRobo_FYIProgram_FunctionNet](#)
- class [CRobo_FYITemp_FunctionNet](#)
- class [CRoboDacqNet](#)
- class [CRoboDeviceNet](#)
CRoboDeviceNet is the base class for all Robo platform based devices
- class [CRoboFluidDeviceNet](#)
- class [CRoboInjectDeviceNet](#)
CRoboInjectDeviceNet is the to control the MCS RoboInject device
- class [CRoboocyte2DeviceNet](#)
CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device
- class [CRoboStatorDeviceNet](#)
- class [CSafeISDeviceNet](#)
- class [CSCUDacqGroupChannelSelectionNet](#)
- class [CSCUFunctionNet](#)
CSCUFunctionNet is the class to control the SCU device
- class [CSerialPortNet](#)
- class [CStg200xBasicNet](#)

- Base class for the Stg200x.*
- class [CStg200xDownloadBasicNet](#)
 - CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.*
- class [CStg200xDownloadNet](#)
 - Main class for the STG download mode This class implements the STG download mode interface.*
- class [CStimulusFunctionNet](#)
- class [CSw2to64DeviceNet](#)
 - The class to control the MCS-USB-Sw2to64 device.*
- class [CTcxDeviceNet](#)
 - Class to control a Temperature Controller (TCX)*
- class [CTEERFunctionNet](#)
 - CTEERFunctionNet is the class to control the TEER device*
- class [CTEERMachineDeviceNet](#)
- class [CUsbExceptionNet](#)
 - Exception class that is thrown in case of an USB error.*
- class [CW2100_FunctionNet](#)
- class [CW2100_StimulatorFunctionNet](#)
- class [CW2100DacqGroupChannelSelectionNet](#)
- class [CWarnerUssingDeviceNet](#)
 - CWarnerUssingDeviceNet is the class to control the Ussing device*
- class [CWarnerUssingFunctionNet](#)
 - CWarnerUssingFunctionNet is the class to control the Ussing device*
- class [CWarnerValveControllerDeviceNet](#)
 - CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller*
- class [CWarnerValveControllerDeviceTesterFunctionNet](#)
 - CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester*
- class [CWClassicFunctionNet](#)
- class [CWirelessBaseFunctionNet](#)
- struct [DeviceldNet](#)
 - Device Id.*
- class [DigitalSource](#)
- class [DriverVersionNet](#)
 - Class gives firmware versions of the device's firmware destinations.*
- class [FirmwareDestinationNames](#)
- class [HeadStageIDType](#)
- class [HeadstageIDTypeObject](#)
- class [HeadStageIDTypeState](#)
- class [mkfilterNet](#)
- class [StgStatusNet](#)
- class [usbSetupPacket_t](#)
- struct [W2100_StimulusParametersNet](#)

Enumerations

- enum [enCMosMeaChipType](#) {
 - [unknown](#) = 0,
 - [nMos16LV](#) = 1,
 - [nMos32LV](#) = 3,
 - [nMos36LN](#) = 6,
 - [nMos64LN](#) = 7 }
- enum [EnSTG200x_STATUS](#) {
 - [OK](#),
 - [NOT_CONNECTED](#),
 - [DEVICE_NOT_FOUND](#) }

Functions

- public delegate void [OnMcsUsbDeviceState](#) (usbSetupPacket_t^ request)
- private delegate void [OnMcsUsbDeviceStateCallback](#) (IntPtr pThis, uint32_t size, IntPtr buffer)
- public delegate void [OnUpdateFirmwareStatusChange](#) (String^)
- public delegate void [OnUpdateFirmwareProgress](#) (int)
- public delegate void [OnDeviceArrivalRemoval](#) (CMcsUsbListEntryNet^ entry)
Delegate to show a device arrival or removal.
- public delegate void [OnStgPollStatus](#) (unsigned int status, [StgStatusNet](#)^ stgStatusNet, array< int >^ index_list)
- public delegate void [OnMwPollStatus](#) (unsigned int CurrentTemp, unsigned int PlateState, unsigned int SwitchState)
- public delegate void [RoboStatusEventDelegate](#) (array< unsigned char >^ buffer)
- public delegate void [OnStg200xDataHandler](#) (uint32_t trigger)
- public delegate void [OnStg200xErrorHandler](#) ()
- public delegate void [OnChannelData](#) (CMcsUsbDacqNet^ dacq, int CbHandle, int numFrames)
- public delegate void [OnError](#) (String^ msg, int action)

9.2.1 Enumeration Type Documentation

9.2.1.1 enCMosMeaChipType `enum enCMosMeaChipType [strong]`

Enumerator

unknown	
nMos16LV	
nMos32LV	
nMos36LN	
nMos64LN	

9.2.1.2 EnSTG200x_STATUS `enum EnSTG200x_STATUS [strong]`

Enumerator

OK	
NOT_CONNECTED	
DEVICE_NOT_FOUND	

9.2.2 Function Documentation

9.2.2.1 OnChannelData() public delegate void Mcs::Usb::OnChannelData (
 CMcsUsbDacqNet[^] *dacq*,
 int *CbHandle*,
 int *numFrames*)

9.2.2.2 OnDeviceArrivalRemoval() public delegate void Mcs::Usb::OnDeviceArrivalRemoval (
 CMcsUsbListEntryNet[^] *entry*)

Delegate to show a device arrival or removal.

9.2.2.3 OnError() public delegate void Mcs::Usb::OnError (
 String[^] *msg*,
 int *action*)

9.2.2.4 OnMcsUsbDeviceState() public delegate void OnMcsUsbDeviceState (
 usbSetupPacket_t[^] *request*)

9.2.2.5 OnMcsUsbDeviceStateCallback() private delegate void OnMcsUsbDeviceStateCallback (
 IntPtr *pThis*,
 uint32_t *size*,
 IntPtr *buffer*)

9.2.2.6 OnMwPollStatus() public delegate void Mcs::Usb::OnMwPollStatus (
 unsigned int *CurrentTemp*,
 unsigned int *PlateState*,
 unsigned int *SwitchState*)

9.2.2.7 OnStg200xDataHandler() public delegate void Mcs::Usb::OnStg200xDataHandler (
 uint32_t *trigger*)

9.2.2.8 OnStg200xErrorHandler() public delegate void Mcs::Usb::OnStg200xErrorHandler ()

9.2.2.9 OnStgPollStatus() `public delegate void Mcs::Usb::OnStgPollStatus (`
`unsigned int status,`
`StgStatusNet^ stgStatusNet,`
`array< int >^ index_list)`

9.2.2.10 OnUpdateFirmwareProgress() `public delegate void Mcs::Usb::OnUpdateFirmwareProgress (`
`int)`

9.2.2.11 OnUpdateFirmwareStatusChange() `public delegate void Mcs::Usb::OnUpdateFirmware↵`
`StatusChange (`
`String^)`

9.2.2.12 RoboStatusEventDelegate() `public delegate void Mcs::Usb::RoboStatusEventDelegate (`
`array< unsigned char >^ buffer)`

10 Class Documentation

10.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference

Public Attributes

- W2100DacqGroupChannelEnumNet [dacqgroup](#)
- int [channel](#)
- int [amplification](#)

10.1.1 Member Data Documentation

10.1.1.1 amplification `int amplification`

10.1.1.2 channel `int channel`

10.1.1.3 dacqgroup `W2100DacqGroupChannelEnumNet dacqgroup`

10.2 BatteryState Class Reference

Properties

- unsigned int [Charge](#) [get]
- unsigned int [Voltage](#) [get]
- System::String^ [ChargeString](#) [get]
- System::String^ [ChargeRegionString](#) [get]
- System::String^ [VoltageString](#) [get]

10.2.1 Property Documentation

10.2.1.1 Charge unsigned int Charge [get]

10.2.1.2 ChargeRegionString System:: String^ ChargeRegionString [get]


10.2.1.3 ChargeString System:: String^ ChargeString [get]

10.2.1.4 Voltage unsigned int Voltage [get]

10.2.1.5 VoltageString System:: String^ VoltageString [get]

10.3 BesselFilterHighPassNet Class Reference

Inheritance diagram for BesselFilterHighPassNet:

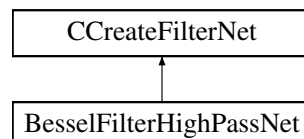

Public Member Functions

- [BesselFilterHighPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)


Additional Inherited Members

10.3.1 Constructor & Destructor Documentation

10.3.1.1 BesselFilterHighPassNet() `BesselFilterHighPassNet (`
 `int numCoefSets,`
 `int order,`
 `double sampleRate,`
 `double cutoffFrequency,`
 `double scale)`

10.4 BesselFilterLowPassNet Class Reference

Inheritance diagram for BesselFilterLowPassNet:

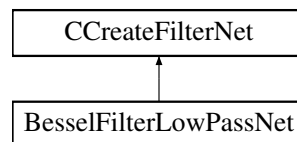

Public Member Functions

- `BesselFilterLowPassNet` (`int numCoefSets`, `int order`, `double sampleRate`, `double cutoffFrequency`, `double scale`)


Additional Inherited Members

10.4.1 Constructor & Destructor Documentation

10.4.1.1 BesselFilterLowPassNet() `BesselFilterLowPassNet (`
 `int numCoefSets,`
 `int order,`
 `double sampleRate,`
 `double cutoffFrequency,`
 `double scale)`

10.5 ButterworthFilterHighPassNet Class Reference

Inheritance diagram for ButterworthFilterHighPassNet:

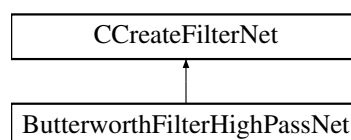

Public Member Functions

- [ButterworthFilterHighPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members

10.5.1 Constructor & Destructor Documentation

10.5.1.1 ButterworthFilterHighPassNet() [ButterworthFilterHighPassNet](#) (


```

    int numCoefSets,
    int order,
    double sampleRate,
    double cutoffFrequency,
    double scale )

```

10.6 ButterworthFilterLowPassNet Class Reference

Inheritance diagram for ButterworthFilterLowPassNet:

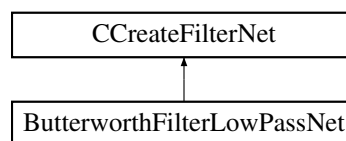

Public Member Functions

- [ButterworthFilterLowPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members

10.6.1 Constructor & Destructor Documentation

10.6.1.1 ButterworthFilterLowPassNet() [ButterworthFilterLowPassNet](#) (


```

    int numCoefSets,
    int order,
    double sampleRate,
    double cutoffFrequency,
    double scale )

```

10.7 CChannelTestDeviceNet Class Reference

Inheritance diagram for CChannelTestDeviceNet:

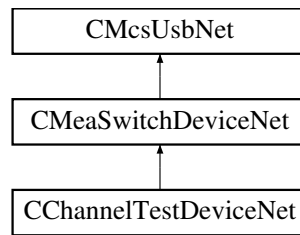

Public Member Functions

- [CChannelTestDeviceNet](#) ()
- [~CChannelTestDeviceNet](#) ()
- void [SetWaveform](#) (unsigned int Waveform)
- void [SetAmplitude](#) (unsigned int Amplitude)
- void [SetFrequency](#) (unsigned int Frequency)
- void [SetAttenuation](#) (unsigned int Attenuation)

Additional Inherited Members

10.7.1 Constructor & Destructor Documentation

10.7.1.1 CChannelTestDeviceNet() [CChannelTestDeviceNet](#) ()

10.7.1.2 ~CChannelTestDeviceNet() [~CChannelTestDeviceNet](#) ()

10.7.2 Member Function Documentation

10.7.2.1 SetAmplitude() void SetAmplitude (
 unsigned int *Amplitude*)


10.7.2.2 SetAttenuation() void SetAttenuation (
 unsigned int *Attenuation*)

10.7.2.3 SetFrequency() `void SetFrequency (`
`unsigned int Frequency)`

10.7.2.4 SetWaveform() `void SetWaveform (`
`unsigned int Waveform)`

10.8 CCMOSMea_FunctionNet Class Reference

Inheritance diagram for CCMOSMea_FunctionNet:

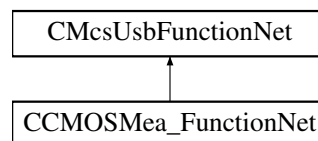

Public Member Functions

- [CCMOSMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] cMOSMea_↔
FunctionPointerContainer)
- [CCMOSMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetADCInputOffset](#) (int32_t offset)
- int32_t [GetADCInputOffset](#) ()
- void [SetSourceDrain](#) (int32_t voltage)
- int32_t [GetSourceDrain](#) ()
- void [SetSourceGate](#) (int32_t voltage)
- int32_t [GetSourceGate](#) ()
- void [SetSourceBulk](#) (int32_t voltage)
- int32_t [GetSourceBulk](#) ()
- void [SetGate](#) (int32_t voltage)
- int32_t [GetGate](#) ()
- void [SetBath](#) (int32_t voltage)
- int32_t [GetBath](#) ()
- int32_t [GetGNDI](#) ()
- int32_t [GetVDDI](#) ()
- int32_t [GetVDD3I](#) ()
- void [UpdateTransistorVoltages](#) ()
- bool [AreTransistorVoltagesSet](#) ()
- void [PowerChip](#) (bool on)
- bool [IsChipPowered](#) ()
- [enCMosMeaChipType](#) [DetectChipType](#) ()
- void [SetGateToVOP](#) ()
- void [SetGateFloating](#) ()
- bool [IsGateFloating](#) ()
- void [VOPSTimerSetResetTimes](#) (uint32_t ResetTime, uint32_t IntervalTime)
- void [VOPSTimerSetResetTimes](#) (uint32_t ResetTime, uint32_t IntervalTime, uint32_t HPFilterResetTime)
- void [SetBathMode](#) (CMOSMeaBathModeEnumNet Mode)
- CMOSMeaBathModeEnumNet [GetBathMode](#) ()
- void [SetNeurochipMemoryData](#) (uint16_t MemAddress, uint32_t MemData)

- void [SetNeurochipMemoryData](#) (uint16_t MemAddress, array< uint32_t >^ MemData)
- uint32_t [GetNeurochipMemoryData](#) (uint16_t MemAddress)
- array< uint32_t > ^ [GetNeurochipMemoryData](#) (uint16_t MemAddress, uint32_t RequestLength)
- uint32_t [GetNeurochipMemorySize](#) ()
- uint32_t [GetMaxNumOfColumns](#) (uint32_t Samplerate)
- void [SetStimulusSites](#) (List< int16_t >^ SwitchPosition)
- List< int16_t > ^ [GetStimulusSites](#) ()
- void [ClearSTGOutput](#) (uint32_t Channel)
- uint32_t [GetNumberOfSupportedGroups](#) ()
- uint32_t [GetNumberOfSupportedGroups](#) (uint32_t virtualDevice)
- DacqGroupChannelEnumNet [GetGroupID](#) (uint32_t Index)
- DacqGroupChannelEnumNet [GetGroupID](#) (uint32_t Index, uint32_t virtualDevice)
- uint32_t [GetGroupNumberOfChannels](#) (DacqGroupChannelEnumNet GroupID)
- uint32_t [GetGroupNumberOfChannels](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- DacqMeaGroupTypeEnumNet [GetGroupType](#) (DacqGroupChannelEnumNet GroupID)
- DacqMeaGroupTypeEnumNet [GetGroupType](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- void [EnableChannelsInGroup](#) (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBitMap)
- void [EnableChannelsInGroup](#) (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBitMap, uint32_t virtualDevice)
- List< bool > ^ [GetEnabledChannelsInGroup](#) (DacqGroupChannelEnumNet GroupID)
- List< bool > ^ [GetEnabledChannelsInGroup](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- SampleSizeNet [GetGroupSampleSize](#) (DacqGroupChannelEnumNet GroupID)
- SampleSizeNet [GetGroupSampleSize](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- uint32_t [GetGroupResolutionPerDigit](#) (DacqGroupChannelEnumNet GroupID)
- uint32_t [GetGroupResolutionPerDigit](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- CMOSMeaValueUnitEnumNet [GetGroupUnit](#) (DacqGroupChannelEnumNet GroupID)
- CMOSMeaValueUnitEnumNet [GetGroupUnit](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- int32_t [GetGroupDCOffset](#) (DacqGroupChannelEnumNet GroupID)
- int32_t [GetGroupDCOffset](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- int32_t [GetGroupADCBits](#) (DacqGroupChannelEnumNet GroupID)
- int32_t [GetGroupADCBits](#) (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- uint32_t [GetGroupChannelBitmaskBySelect](#) (DacqGroupChannelEnumNet GroupID, uint32_t ChannelNumber)
- uint32_t [GetGroupChannelBitmaskBySelect](#) (DacqGroupChannelEnumNet GroupID, uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaInterfaceADCEnumNet [GetGroupChannelBitmaskInterfaceADC](#) (uint32_t ChannelNumber)
- CMOSMeaInterfaceADCEnumNet [GetGroupChannelBitmaskInterfaceADC](#) (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaIFDigChannelEnumNet [GetGroupChannelBitmaskIFDigChannels](#) (uint32_t ChannelNumber)
- CMOSMeaIFDigChannelEnumNet [GetGroupChannelBitmaskIFDigChannels](#) (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCBathCurrentEnumNet [GetGroupChannelBitmaskHS1NCBathCurrent](#) (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCBathCurrentEnumNet [GetGroupChannelBitmaskHS1NCBathCurrent](#) (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet [GetGroupChannelBitmaskHS1NCCol2Current](#) (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet [GetGroupChannelBitmaskHS1NCCol2Current](#) (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NChipTempEnumNet [GetGroupChannelBitmaskHS1NChipTemp](#) (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NChipTempEnumNet [GetGroupChannelBitmaskHS1NChipTemp](#) (uint32_t ChannelNumber, uint32_t virtualDevice)

- CMOSMeaSTG1DACSignalEnumNet [GetGroupChannelBitmaskSTG1DACSignal](#) (uint32_t Channel↔ Number)
- CMOSMeaSTG1DACSignalEnumNet [GetGroupChannelBitmaskSTG1DACSignal](#) (uint32_t Channel↔ Number, uint32_t virtualDevice)
- CMOSMeaHS1SidebandEnumNet [GetGroupChannelBitmaskHS1Sidebands](#) (uint32_t ChannelNumber)
- CMOSMeaHS1SidebandEnumNet [GetGroupChannelBitmaskHS1Sidebands](#) (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHS1TriggerStatusEnumNet [GetGroupChannelBitmaskHS1TriggerStatus](#) (uint32_t Channel↔ Number)
- CMOSMeaHS1TriggerStatusEnumNet [GetGroupChannelBitmaskHS1TriggerStatus](#) (uint32_t Channel↔ Number, uint32_t virtualDevice)
- CMOSMeaPacketFrameContextGroupEnumNet [GetGroupChannelBitmaskPacketFrameContext](#) (uint32_t↔ ChannelNumber)
- CMOSMeaPacketFrameContextGroupEnumNet [GetGroupChannelBitmaskPacketFrameContext](#) (uint32_t↔ ChannelNumber, uint32_t virtualDevice)

Additional Inherited Members

10.8.1 Constructor & Destructor Documentation

10.8.1.1 CCMOSMea_FunctionNet() [1/2] [CCMOSMea_FunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *cCMOSMea_FunctionPointerContainer*)

10.8.1.2 CCMOSMea_FunctionNet() [2/2] [CCMOSMea_FunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

10.8.2 Member Function Documentation

10.8.2.1 AreTransistorVoltagesSet() [bool AreTransistorVoltagesSet](#) ()

10.8.2.2 ClearSTGOutput() [void ClearSTGOutput](#) (
 uint32_t *Channel*)

10.8.2.3 DetectChipType() [enCMosMeaChipType DetectChipType](#) ()

10.8.2.4 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (
DacqGroupChannelEnumNet GroupID,
List< bool >^ EnabledChannelsBitMap)

10.8.2.5 EnableChannelsInGroup() [2/2] void EnableChannelsInGroup (
DacqGroupChannelEnumNet GroupID,
List< bool >^ EnabledChannelsBitMap,
uint32_t virtualDevice)

10.8.2.6 GetADCInputOffset() int32_t GetADCInputOffset ()

10.8.2.7 GetBath() int32_t GetBath ()

10.8.2.8 GetBathMode() CMOSMeaBathModeEnumNet GetBathMode ()

10.8.2.9 GetEnabledChannelsInGroup() [1/2] List<bool> ^ GetEnabledChannelsInGroup (
DacqGroupChannelEnumNet GroupID)

10.8.2.10 GetEnabledChannelsInGroup() [2/2] List<bool> ^ GetEnabledChannelsInGroup (
DacqGroupChannelEnumNet GroupID,
uint32_t virtualDevice)

10.8.2.11 GetGate() int32_t GetGate ()

10.8.2.12 GetGNDI() int32_t GetGNDI ()

10.8.2.13 GetGroupADCBits() [1/2] int32_t GetGroupADCBits (
DacqGroupChannelEnumNet GroupID)

10.8.2.14 GetGroupADCBits() [2/2] `int32_t GetGroupADCBits (`
`DacqGroupChannelEnumNet GroupID,`
`uint32_t virtualDevice)`

10.8.2.15 GetGroupChannelBitmaskBySelect() [1/2] `uint32_t GetGroupChannelBitmaskBySelect (`
`DacqGroupChannelEnumNet GroupID,`
`uint32_t ChannelNumber)`

10.8.2.16 GetGroupChannelBitmaskBySelect() [2/2] `uint32_t GetGroupChannelBitmaskBySelect (`
`DacqGroupChannelEnumNet GroupID,`
`uint32_t ChannelNumber,`
`uint32_t virtualDevice)`

10.8.2.17 GetGroupChannelBitmaskHS1NCBathCurrent() [1/2] `CMOSMeaHeadstage1NCBathCurrentEnum↔`
`Net GetGroupChannelBitmaskHS1NCBathCurrent (`
`uint32_t ChannelNumber)`

10.8.2.18 GetGroupChannelBitmaskHS1NCBathCurrent() [2/2] `CMOSMeaHeadstage1NCBathCurrentEnum↔`
`Net GetGroupChannelBitmaskHS1NCBathCurrent (`
`uint32_t ChannelNumber,`
`uint32_t virtualDevice)`

10.8.2.19 GetGroupChannelBitmaskHS1NCCol2Current() [1/2] `CMOSMeaHeadstage1NCCol2CurrentEnum↔`
`Net GetGroupChannelBitmaskHS1NCCol2Current (`
`uint32_t ChannelNumber)`

10.8.2.20 GetGroupChannelBitmaskHS1NCCol2Current() [2/2] `CMOSMeaHeadstage1NCCol2CurrentEnum↔`
`Net GetGroupChannelBitmaskHS1NCCol2Current (`
`uint32_t ChannelNumber,`
`uint32_t virtualDevice)`

10.8.2.21 GetGroupChannelBitmaskHS1NChipTemp() [1/2] `CMOSMeaHeadstage1NChipTempEnumNet Get↔`
`GroupChannelBitmaskHS1NChipTemp (`
`uint32_t ChannelNumber)`

10.8.2.22 GetGroupChannelBitmaskHS1NChipTemp() [2/2] CMOSMeaHeadstage1NChipTempEnumNet Get↔
 GroupChannelBitmaskHS1NChipTemp (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

10.8.2.23 GetGroupChannelBitmaskHS1Sidebands() [1/2] CMOSMeaHS1SidebandEnumNet GetGroup↔
 ChannelBitmaskHS1Sidebands (
 uint32_t ChannelNumber)

10.8.2.24 GetGroupChannelBitmaskHS1Sidebands() [2/2] CMOSMeaHS1SidebandEnumNet GetGroup↔
 ChannelBitmaskHS1Sidebands (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

10.8.2.25 GetGroupChannelBitmaskHS1TriggerStatus() [1/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
 GroupChannelBitmaskHS1TriggerStatus (
 uint32_t ChannelNumber)

10.8.2.26 GetGroupChannelBitmaskHS1TriggerStatus() [2/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
 GroupChannelBitmaskHS1TriggerStatus (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

10.8.2.27 GetGroupChannelBitmaskIFDigChannels() [1/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
 ChannelBitmaskIFDigChannels (
 uint32_t ChannelNumber)

10.8.2.28 GetGroupChannelBitmaskIFDigChannels() [2/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
 ChannelBitmaskIFDigChannels (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

10.8.2.29 GetGroupChannelBitmaskInterfaceADC() [1/2] CMOSMeaInterfaceADCEnumNet GetGroup↔
 ChannelBitmaskInterfaceADC (
 uint32_t ChannelNumber)

10.8.2.30 GetGroupChannelBitmaskInterfaceADC() [2/2] CMOSMeaInterfaceADCEnumNet GetGroup↔
 ChannelBitmaskInterfaceADC (
 uint32_t *ChannelNumber*,
 uint32_t *virtualDevice*)

10.8.2.31 GetGroupChannelBitmaskPacketFrameContext() [1/2] CMOSMeaPacketFrameContextGroup↔
 EnumNet GetGroupChannelBitmaskPacketFrameContext (
 uint32_t *ChannelNumber*)

10.8.2.32 GetGroupChannelBitmaskPacketFrameContext() [2/2] CMOSMeaPacketFrameContextGroup↔
 EnumNet GetGroupChannelBitmaskPacketFrameContext (
 uint32_t *ChannelNumber*,
 uint32_t *virtualDevice*)

10.8.2.33 GetGroupChannelBitmaskSTG1DACSignal() [1/2] CMOSMeaSTG1DACSignalEnumNet GetGroup↔
 ChannelBitmaskSTG1DACSignal (
 uint32_t *ChannelNumber*)

10.8.2.34 GetGroupChannelBitmaskSTG1DACSignal() [2/2] CMOSMeaSTG1DACSignalEnumNet GetGroup↔
 ChannelBitmaskSTG1DACSignal (
 uint32_t *ChannelNumber*,
 uint32_t *virtualDevice*)

10.8.2.35 GetGroupDCOffset() [1/2] int32_t GetGroupDCOffset (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.36 GetGroupDCOffset() [2/2] int32_t GetGroupDCOffset (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.37 GetGroupID() [1/2] DacqGroupChannelEnumNet GetGroupID (
 uint32_t *Index*)

10.8.2.38 GetGroupID() [2/2] DacqGroupChannelEnumNet GetGroupID (
 uint32_t *Index*,
 uint32_t *virtualDevice*)

10.8.2.39 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.40 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.41 GetGroupResolutionPerDigit() [1/2] uint32_t GetGroupResolutionPerDigit (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.42 GetGroupResolutionPerDigit() [2/2] uint32_t GetGroupResolutionPerDigit (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.43 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.44 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.45 GetGroupType() [1/2] DacqMeaGroupTypeEnumNet GetGroupType (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.46 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.47 GetGroupUnit() [1/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
 DacqGroupChannelEnumNet *GroupID*)

10.8.2.48 GetGroupUnit() [2/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
 DacqGroupChannelEnumNet *GroupID*,
 uint32_t *virtualDevice*)

10.8.2.49 GetMaxNumOfColumns() uint32_t GetMaxNumOfColumns (
 uint32_t *Samplerate*)

10.8.2.50 GetNeurochipMemoryData() [1/2] uint32_t GetNeurochipMemoryData (
 uint16_t *MemAddress*)

10.8.2.51 GetNeurochipMemoryData() [2/2] array<uint32_t> ^ GetNeurochipMemoryData (
 uint16_t *MemAddress*,
 uint32_t *RegestLength*)

10.8.2.52 GetNeurochipMemorySize() uint32_t GetNeurochipMemorySize ()

10.8.2.53 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ()

10.8.2.54 GetNumberOfSupportedGroups() [2/2] uint32_t GetNumberOfSupportedGroups (
 uint32_t *virtualDevice*)

10.8.2.55 GetSourceBulk() int32_t GetSourceBulk ()

10.8.2.56 GetSourceDrain() int32_t GetSourceDrain ()

10.8.2.57 GetSourceGate() `int32_t GetSourceGate ()`

10.8.2.58 GetStimulusSites() `List<int16_t> ^ GetStimulusSites ()`

10.8.2.59 GetVDD3I() `int32_t GetVDD3I ()`

10.8.2.60 GetVDDI() `int32_t GetVDDI ()`

10.8.2.61 IsChipPowered() `bool IsChipPowered ()`

10.8.2.62 IsGateFloating() `bool IsGateFloating ()`

10.8.2.63 PowerChip() `void PowerChip (`
`bool on)`

10.8.2.64 SetADCInputOffset() `void SetADCInputOffset (`
`int32_t offset)`

10.8.2.65 SetBath() `void SetBath (`
`int32_t voltage)`

10.8.2.66 SetBathMode() `void SetBathMode (`
`CMOSMeaBathModeEnumNet Mode)`

10.8.2.67 SetGate() void SetGate (
 int32_t voltage)

10.8.2.68 SetGateFloating() void SetGateFloating ()

10.8.2.69 SetGateToVOP() void SetGateToVOP ()

10.8.2.70 SetNeurochipMemoryData() [1/2] void SetNeurochipMemoryData (
 uint16_t MemAddress,
 array< uint32_t >^ MemData)

10.8.2.71 SetNeurochipMemoryData() [2/2] void SetNeurochipMemoryData (
 uint16_t MemAddress,
 uint32_t MemData)

10.8.2.72 SetSourceBulk() void SetSourceBulk (
 int32_t voltage)

10.8.2.73 SetSourceDrain() void SetSourceDrain (
 int32_t voltage)

10.8.2.74 SetSourceGate() void SetSourceGate (
 int32_t voltage)

10.8.2.75 SetStimulusSites() void SetStimulusSites (
 List< int16_t >^ SwitchPosition)


10.8.2.76 UpdateTransistorVoltages() void UpdateTransistorVoltages ()

10.8.2.77 VOPSTimerSetResetTimes() [1/2] `void VOPSTimerSetResetTimes (`
`uint32_t ResetTime,`
`uint32_t IntervalTime)`

10.8.2.78 VOPSTimerSetResetTimes() [2/2] `void VOPSTimerSetResetTimes (`
`uint32_t ResetTime,`
`uint32_t IntervalTime,`
`uint32_t HPFilterResetTime)`

10.9 CCMOSMeaDeviceNet Class Reference

Inheritance diagram for CCMOSMeaDeviceNet:

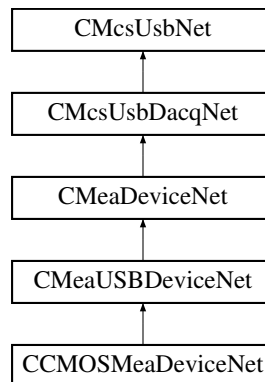

Classes

- class [CRegionOfInterestRect](#)

Public Member Functions

- [CCMOSMeaDeviceNet](#) (void)
- [~CCMOSMeaDeviceNet](#) ()
- void [SetBaseSamplerate](#) (int BaseSamplerate)
- int [GetBaseSamplerate](#) ()
- array< int > ^ [GetAvailableBaseSamplerates](#) ()
- void [SetRegionOfInterests](#) (System::Collections::Generic::Dictionary< int, [CRegionOfInterestRect](#)>^ rois)
- void [UpdateChannelBlock](#) (int queuesize, int threshold, int channels_in_block)
- System::Collections::Generic::Dictionary< int, array< array< int16_t >^>^ > ^ [GetCMOSDataDictionary](#) (int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint16_t >^>^ > ^ [GetChannelDataUI16](#) (DacqGroup< ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int16_t >^>^ > ^ [GetChannelDataI16](#) (DacqGroup< ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint32_t >^>^ > ^ [GetChannelDataUI32](#) (DacqGroup< ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int32_t >^>^ > ^ [GetChannelDataI32](#) (DacqGroup< ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Properties

- [CCMOSMea_FunctionNet](#)^ [CMosMea](#) [get]
- [CStimulusFunctionNet](#)^ [Stimulus](#) [get]

Additional Inherited Members

10.9.1 Constructor & Destructor Documentation

10.9.1.1 CCMOSMeaDeviceNet() [CCMOSMeaDeviceNet](#) (
 void)

10.9.1.2 ~CCMOSMeaDeviceNet() [~CCMOSMeaDeviceNet](#) ()

10.9.2 Member Function Documentation

10.9.2.1 GetAvailableBaseSamplerates() [array<int>](#) ^ [GetAvailableBaseSamplerates](#) ()

10.9.2.2 GetBaseSamplerate() [int](#) [GetBaseSamplerate](#) ()

10.9.2.3 GetChannelDataI16() [System::Collections::Generic::Dictionary<int, array<int16_t>^>](#) ^ [GetChannelDataI16](#) (
 [DacqGroupChannelEnumNet](#) *group*,
 [int](#) *frames*,
 [System::Runtime::InteropServices::Out] [int](#) % *frames_ret*)

10.9.2.4 GetChannelDataI32() [System::Collections::Generic::Dictionary<int, array<int32_t>^>](#) ^ [GetChannelDataI32](#) (
 [DacqGroupChannelEnumNet](#) *group*,
 [int](#) *frames*,
 [System::Runtime::InteropServices::Out] [int](#) % *frames_ret*)

10.9.2.5 GetChannelDataUI16() `System::Collections::Generic::Dictionary<int, array<uint16_t>^>`
`^ GetChannelDataUI16 (`
`DacqGroupChannelEnumNet group,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret)`

10.9.2.6 GetChannelDataUI32() `System::Collections::Generic::Dictionary<int, array<uint32_t>^>`
`^ GetChannelDataUI32 (`
`DacqGroupChannelEnumNet group,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret)`

10.9.2.7 GetCMOSDataDictionary() `System::Collections::Generic::Dictionary<int, array<array<int16_t>^>^>`
`^ GetCMOSDataDictionary (`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret)`

10.9.2.8 SetBaseSamplerate() `void SetBaseSamplerate (`
`int BaseSamplerate)`

10.9.2.9 SetRegionOfInterests() `void SetRegionOfInterests (`
`System::Collections::Generic::Dictionary< int, CRegionOfInterestRect^>^ rois)`

10.9.2.10 UpdateChannelBlock() `void UpdateChannelBlock (`
`int queuesize,`
`int threshold,`
`int channels_in_block)`


10.9.3 Property Documentation

10.9.3.1 CMosMea `CCMOSMea_FunctionNet^ CMosMea [get]`

10.9.3.2 Stimulus `CStimulusFunctionNet^ Stimulus [get]`

10.10 CCreateFilterNet Class Reference

Inheritance diagram for CCreateFilterNet:

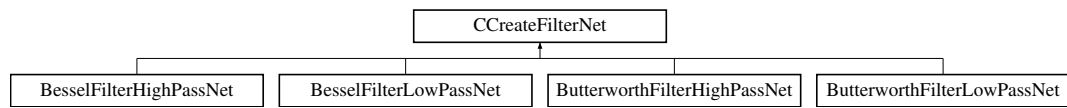

Public Member Functions

- [CCreateFilterNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)
- [~CCreateFilterNet](#) ()
- [CFilterCoefficientsNet](#) ^ [GetBiQuad](#) (int index)
- array< [CFilterCoefficientsNet](#) ^> ^ [GetBiQuads](#) ()

Static Public Member Functions

- static int [FindFilter](#) (array< [CFilterCoefficientsNet](#) ^>^ coef, array< [CCreateFilterNet](#) ^>^ param)

Protected Member Functions

- [CCreateFilterNet](#) (int numCoefSets, CCreateFilter *pCreateFilter)

Properties

- int [NumCoefSets](#) [get]
- int [Order](#) [get]
- double [SampleRate](#) [get]
- double [CutoffFrequency](#) [get]
- double [Scale](#) [get]

10.10.1 Constructor & Destructor Documentation

10.10.1.1 CCreateFilterNet() [1/2] [CCreateFilterNet](#) (

```

    int numCoefSets,
    int order,
    double sampleRate,
    double cutoffFrequency,
    double scale )
  
```

10.10.1.2 ~CCreateFilterNet() [~CCreateFilterNet](#) ()

10.10.1.3 CCreateFilterNet() [2/2] `CCreateFilterNet` (
 int numCoefSets,
 `CCreateFilter` * *pCreateFilter*) [protected]

10.10.2 Member Function Documentation

10.10.2.1 FindFilter() `static int FindFilter` (
 array< `CFilterCoefficientsNet` ^>^ *coef*,
 array< `CCreateFilterNet` ^>^ *param*) [static]

10.10.2.2 GetBiQuad() `CFilterCoefficientsNet` ^ GetBiQuad (
 int index)

10.10.2.3 GetBiQuads() array<`CFilterCoefficientsNet` ^> ^ GetBiQuads ()

10.10.3 Property Documentation

10.10.3.1 CutoffFrequency `double CutoffFrequency` [get]

10.10.3.2 NumCoefSets `int NumCoefSets` [get]


10.10.3.3 Order `int Order` [get]

10.10.3.4 SampleRate `double SampleRate` [get]

10.10.3.5 Scale `double Scale` [get]

10.11 CDacCalibrationFunctionNet Class Reference

Inheritance diagram for CDacCalibrationFunctionNet:

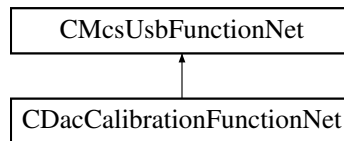

Public Member Functions

- [CDacCalibrationFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pDacCalibrationFunctionPointerContainer)
Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.
- [CDacCalibrationFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CDacCalibrationFunctionNet](#) ()
- [!CDacCalibrationFunctionNet](#) ()
- void [SetDacOffset](#) (uint16_t dacChannel, uint32_t offset)
Sets the offset of a DAC channel.
- uint32_t [GetDacOffset](#) (uint16_t dacChannel)
Gets the offset of a DAC channel.
- void [BurnDacOffset](#) (uint16_t dacChannel)
Writes the offset of a DAC channel to permanent memory.

Additional Inherited Members

10.11.1 Detailed Description

10.11.2 Constructor & Destructor Documentation

10.11.2.1 CDacCalibrationFunctionNet() [1/2] [CDacCalibrationFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pDacCalibrationFunctionPointerContainer)

Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.

10.11.2.2 CDacCalibrationFunctionNet() [2/2] [CDacCalibrationFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.11.2.3 ~CDacCalibrationFunctionNet() virtual [~CDacCalibrationFunctionNet](#) () [virtual]

10.11.2.4 `!CDacCalibrationFunctionNet()` `!CDacCalibrationFunctionNet ()`

10.11.3 Member Function Documentation

10.11.3.1 **BurnDacOffset()** `void BurnDacOffset (` `uint16_t dacChannel)`

Writes the offset of a DAC channel to permanent memory.

Parameters

<i>dacChannel</i>	The DAC channel number.
-------------------	-------------------------

10.11.3.2 **GetDacOffset()** `uint32_t GetDacOffset (` `uint16_t dacChannel)`

Gets the offset of a DAC channel.

Parameters

<i>dacChannel</i>	The DAC channel number.
-------------------	-------------------------

Returns

The offset in digits.

10.11.3.3 **SetDacOffset()** `void SetDacOffset (` `uint16_t dacChannel,` `uint32_t offset)`


Sets the offset of a DAC channel.

Parameters

<i>dacChannel</i>	The DAC channel number.
<i>offset</i>	The offset in digits.

10.12 CDacqGroupChannelGenericSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelGenericSelectionNet:

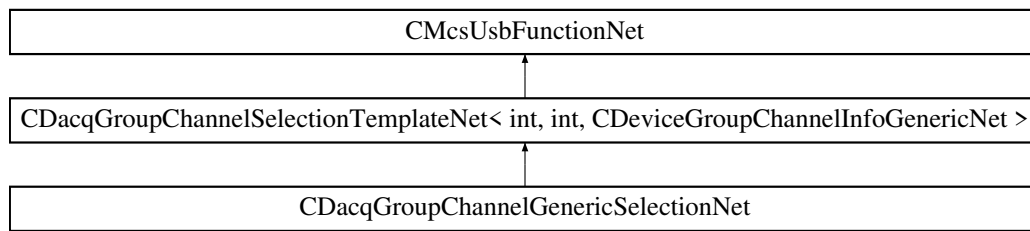

Public Member Functions

- [CDacqGroupChannelGenericSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)


Additional Inherited Members

10.12.1 Constructor & Destructor Documentation

10.12.1.1 CDacqGroupChannelGenericSelectionNet() [CDacqGroupChannelGenericSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)

10.13 CDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelSelectionNet:

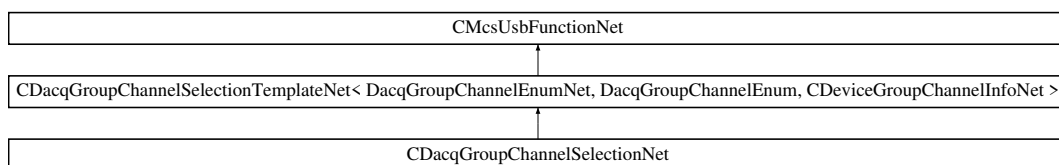

Public Member Functions

- [CDacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)


Additional Inherited Members

10.13.1 Constructor & Destructor Documentation

10.13.1.1 CDacqGroupChannelSelectionNet() [CDacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)

10.14 CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference

Inheritance diagram for CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >:

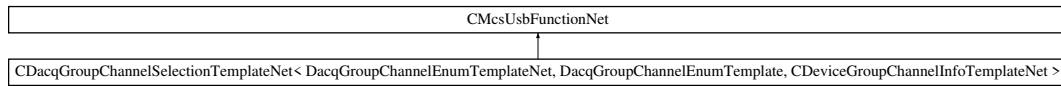

Public Member Functions

- [CDacqGroupChannelSelectionTemplateNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- [uint32_t GetNumberOfSupportedGroups](#) ()
- [uint32_t GetNumberOfSupportedGroups](#) ([uint32_t](#) virtualDevice)
- [DacqGroupChannelEnumTemplateNet GetGroupID](#) ([uint32_t](#) Index)
- [DacqGroupChannelEnumTemplateNet GetGroupID](#) ([uint32_t](#) Index, [uint32_t](#) virtualDevice)
- [uint32_t GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [uint32_t GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [DacqMeaGroupTypeEnumNet GetGroupType](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [DacqMeaGroupTypeEnumNet GetGroupType](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap)
- [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap, [uint32_t](#) virtualDevice)
- [List< bool >](#)[^] [GetEnabledChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [List< bool >](#)[^] [GetEnabledChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [SampleSizeNet GetGroupSampleSize](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [SampleSizeNet GetGroupSampleSize](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [List< CDeviceGroupChannelInfoTemplateNet[^]>](#)[^] [GetDeviceGroupChannelInfos](#) ()
- [List< CDeviceGroupChannelInfoTemplateNet[^]>](#)[^] [GetDeviceGroupChannelInfos](#) ([uint32_t](#) virtualDevice)

Additional Inherited Members

10.14.1 Constructor & Destructor Documentation

10.14.1.1 CDacqGroupChannelSelectionTemplateNet() [CDacqGroupChannelSelectionTemplateNet](#) ([CMcsUsbNet](#)[^] mcsusb)

10.14.2 Member Function Documentation

10.14.2.1 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID,
List< bool >^ EnabledChannelsBitMap)

10.14.2.2 EnableChannelsInGroup() [2/2] void EnableChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID,
List< bool >^ EnabledChannelsBitMap,
uint32_t virtualDevice)

10.14.2.3 GetDeviceGroupChannelInfos() [1/2] List<CDeviceGroupChannelInfoTemplateNet^> ^ Get↔
DeviceGroupChannelInfos ()

10.14.2.4 GetDeviceGroupChannelInfos() [2/2] List<CDeviceGroupChannelInfoTemplateNet^> ^ Get↔
DeviceGroupChannelInfos (
uint32_t virtualDevice)

10.14.2.5 GetEnabledChannelsInGroup() [1/2] List<bool> ^ GetEnabledChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID)

10.14.2.6 GetEnabledChannelsInGroup() [2/2] List<bool> ^ GetEnabledChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID,
uint32_t virtualDevice)

10.14.2.7 GetGroupID() [1/2] DacqGroupChannelEnumTemplateNet GetGroupID (
uint32_t Index)

10.14.2.8 GetGroupID() [2/2] DacqGroupChannelEnumTemplateNet GetGroupID (
uint32_t Index,
uint32_t virtualDevice)

10.14.2.9 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
DacqGroupChannelEnumTemplateNet GroupID)

10.14.2.10 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

10.14.2.11 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
 DacqGroupChannelEnumTemplateNet GroupID)

10.14.2.12 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

10.14.2.13 GetGroupType() [1/2] DacqMeaGroupTypeEnumNet GetGroupType (
 DacqGroupChannelEnumTemplateNet GroupID)


10.14.2.14 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

10.14.2.15 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ()

10.14.2.16 GetNumberOfSupportedGroups() [2/2] uint32_t GetNumberOfSupportedGroups (
 uint32_t virtualDevice)

10.15 CDeviceGroupChannelInfoGenericNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoGenericNet:

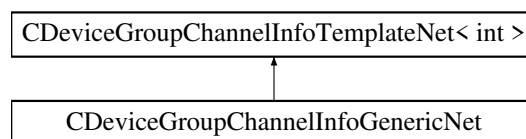

Public Member Functions

- [CDeviceGroupChannelInfoGenericNet](#) (int id, int channels, DacqMeaGroupTypeEnumNet type)


Additional Inherited Members

10.15.1 Constructor & Destructor Documentation

10.15.1.1 CDeviceGroupChannelInfoGenericNet() [CDeviceGroupChannelInfoGenericNet](#) (
 int *id*,
 int *channels*,
 DacqMeaGroupTypeEnumNet *type*)

10.16 CDeviceGroupChannelInfoNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoNet:

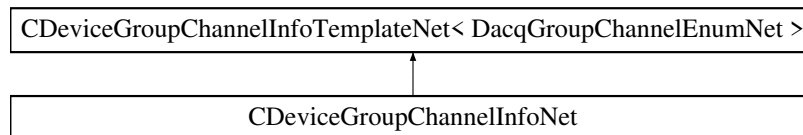

Public Member Functions

- [CDeviceGroupChannelInfoNet](#) (DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)


Additional Inherited Members

10.16.1 Constructor & Destructor Documentation

10.16.1.1 CDeviceGroupChannelInfoNet() [CDeviceGroupChannelInfoNet](#) (
 DacqGroupChannelEnumNet *id*,
 int *channels*,
 DacqMeaGroupTypeEnumNet *type*)

10.17 CDeviceGroupChannelInfoSCUNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoSCUNet:

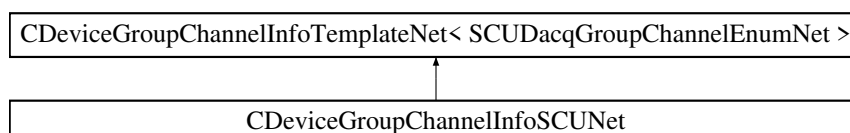

Public Member Functions

- [CDeviceGroupChannelInfoSCUNet](#) (SCUDacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

10.17.1 Constructor & Destructor Documentation

10.17.1.1 CDeviceGroupChannelInfoSCUNet() [CDeviceGroupChannelInfoSCUNet](#) (SCUDacqGroupChannelEnumNet *id*, int *channels*, DacqMeaGroupTypeEnumNet *type*)

10.18 CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet > Class Template Reference

Public Member Functions

- [CDeviceGroupChannelInfoTemplateNet](#) (DacqGroupChannelEnumTemplateNet id, int channels, DacqMeaGroupTypeEnumNet type)

Public Attributes

- DacqGroupChannelEnumTemplateNet [GroupID](#)
- int [NumberOfChannels](#)
- DacqMeaGroupTypeEnumNet [GroupType](#)

10.18.1 Constructor & Destructor Documentation

10.18.1.1 CDeviceGroupChannelInfoTemplateNet() [CDeviceGroupChannelInfoTemplateNet](#) (DacqGroupChannelEnumTemplateNet *id*, int *channels*, DacqMeaGroupTypeEnumNet *type*)

10.18.2 Member Data Documentation


10.18.2.1 GroupID DacqGroupChannelEnumTemplateNet GroupID

10.18.2.2 GroupType `DacqMeaGroupTypeEnumNet GroupType`

10.18.2.3 NumberOfChannels `int NumberOfChannels`

10.19 CDeviceGroupChannelInfoW2100Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoW2100Net:

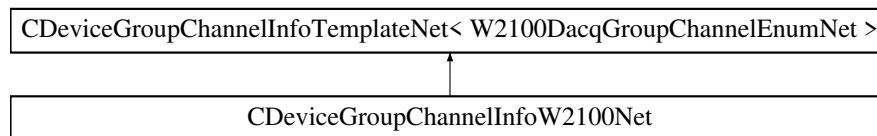

Public Member Functions

- [CDeviceGroupChannelInfoW2100Net](#) ([W2100DacqGroupChannelEnumNet](#) id, int channels, [DacqMeaGroupTypeEnumNet](#) type)

Additional Inherited Members


10.19.1 Constructor & Destructor Documentation

10.19.1.1 CDeviceGroupChannelInfoW2100Net() [CDeviceGroupChannelInfoW2100Net](#) (
 [W2100DacqGroupChannelEnumNet](#) id,
 int channels,
 [DacqMeaGroupTypeEnumNet](#) type)

10.20 CDigOutStimulatorFunctionNet Class Reference

[CDigOutStimulatorFunctionNet](#) is the class of the DigOut stimulator function class.

Inheritance diagram for CDigOutStimulatorFunctionNet:

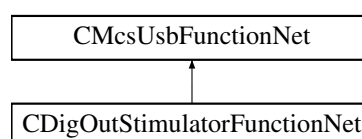

Public Member Functions

- [CDigOutStimulatorFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pDigOutStimulatorFunctionPointerContainer)
Initializes a new instance of the [CDigOutStimulatorFunctionNet](#) class.
- [CDigOutStimulatorFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CDigOutStimulatorFunctionNet](#) ()
- [!CDigOutStimulatorFunctionNet](#) ()
- void [ClearChannel](#) (int32_t NrChannel)
clear stimulation pattern
- [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#) ^ [PrepareChannelData](#) (array< int32_t >^ Amplitude, array< uint64_t >^ Duration)
- void [SendChannelData](#) (int32_t NrChannel, [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] device_data_and_unrolled)
- int32_t [GetNumberOfChannels](#) ()
- void [SetGlobalRepeat](#) (int32_t NrChannel, bool value)
set repeat whole stimulation pattern
- bool [GetGlobalRepeat](#) (int32_t NrChannel)
get repeat whole stimulation pattern
- void [SetStartTriggerSlope](#) (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)
sets start condition of digital out stimulator
- DigitalStimulatorTriggerSlopeEnumNet [GetStartTriggerSlope](#) (int32_t NrChannel)
queries start condition of digital out stimulator
- void [SetStopTriggerSlope](#) (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)
sets stop condition of digital out stimulator
- DigitalStimulatorTriggerSlopeEnumNet [GetStopTriggerSlope](#) (int32_t NrChannel)
queries stop condition of digital out stimulator

Additional Inherited Members**10.20.1 Detailed Description**

[CDigOutStimulatorFunctionNet](#) is the class of the DigOut stimulator function class.

10.20.2 Constructor & Destructor Documentation

- 10.20.2.1 [CDigOutStimulatorFunctionNet](#)() [1/2]** [CDigOutStimulatorFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pDigOutStimulatorFunctionPointerContainer)

Initializes a new instance of the [CDigOutStimulatorFunctionNet](#) class.

- 10.20.2.2 [CDigOutStimulatorFunctionNet](#)() [2/2]** [CDigOutStimulatorFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.20.2.3 `~CDigOutStimulatorFunctionNet()` `virtual ~CDigOutStimulatorFunctionNet () [virtual]`

10.20.2.4 `!CDigOutStimulatorFunctionNet()` `!CDigOutStimulatorFunctionNet ()`

10.20.3 Member Function Documentation

10.20.3.1 `ClearChannel()` `void ClearChannel (`
`int32_t NrChannel)`

clear stimulation pattern

Parameters

<i>NrChannel</i>	
------------------	--

10.20.3.2 `GetGlobalRepeat()` `bool GetGlobalRepeat (`
`int32_t NrChannel)`

get repeat whole stimulation pattern

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

current value

10.20.3.3 `GetNumberOfChannels()` `int32_t GetNumberOfChannels ()`

Returns

10.20.3.4 `GetStartTriggerSlope()` `DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (`
`int32_t NrChannel)`

queries start condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

start condition (rising or falling edge)

10.20.3.5 GetStopTriggerSlope() `DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope (int32_t NrChannel)`

queries stop condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

stop condition (rising or falling edge)

10.20.3.6 PrepareChannelData() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareChannelData (array< int32_t >^ Amplitude, array< uint64_t >^ Duration)`

Parameters

<i>Amplitude</i>	
<i>Duration</i>	

Returns

10.20.3.7 SendChannelData() `void SendChannelData (int32_t NrChannel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled)`

Parameters

<i>NrChannel</i>	
<i>device_data_and_unrolled</i>	

10.20.3.8 SetGlobalRepeat() void SetGlobalRepeat (
 int32_t *NrChannel*,
 bool *value*)

set repeat whole stimulation pattern

Parameters

<i>NrChannel</i>	channel number
<i>value</i>	new value

10.20.3.9 SetStartTriggerSlope() void SetStartTriggerSlope (
 int32_t *NrChannel*,
 DigitalStimulatorTriggerSlopeEnumNet *Condition*)

sets start condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
<i>Condition</i>	start condition (rising or falling edge)

10.20.3.10 SetStopTriggerSlope() void SetStopTriggerSlope (
 int32_t *NrChannel*,
 DigitalStimulatorTriggerSlopeEnumNet *Condition*)

sets stop condition of digital out stimulator


Parameters

<i>NrChannel</i>	channel number
<i>Condition</i>	stop condition (rising or falling edge)

10.21 CEncapsulatorDeviceNet Class Reference

[CEncapsulatorDeviceNet](#) is the to control the MCS HiClamp device

Inheritance diagram for CEncapsulatorDeviceNet:

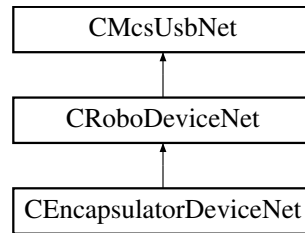

Public Member Functions

- [CEncapsulatorDeviceNet](#) (void)
- [CRoboFluidDeviceNet](#) ^ [GetRoboFluidDevice](#) ()

Additional Inherited Members

10.21.1 Detailed Description

[CEncapsulatorDeviceNet](#) is the to control the MCS HiClamp device

10.21.2 Constructor & Destructor Documentation

10.21.2.1 CEncapsulatorDeviceNet() [CEncapsulatorDeviceNet](#) (
void)


10.21.3 Member Function Documentation

10.21.3.1 GetRoboFluidDevice() [CRoboFluidDeviceNet](#) ^ [GetRoboFluidDevice](#) ()

10.22 CExternDTesterDeviceNet Class Reference

[CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)

Inheritance diagram for CExternDTesterDeviceNet:

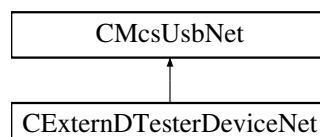

Public Member Functions

- [CExternDTesterDeviceNet](#) ()
Initializes a new instance of the [CExternDTesterDeviceNet](#) class.
- virtual [~CExternDTesterDeviceNet](#) ()
- [IExternDTesterDeviceNet](#) ()
- array< uint8_t > ^ [Read](#) (int configString_Length)
Reads the config string from the device.
- String ^ [Read2](#) ()
Reads the config string from the device.
- void [Write](#) (array< uint8_t > ^ configString)
Reads the config string from the device.
- void [Write2](#) (String ^ configString)
Reads the config string from the device.

Additional Inherited Members**10.22.1 Detailed Description**

[CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)

10.22.2 Constructor & Destructor Documentation**10.22.2.1 CExternDTesterDeviceNet() [CExternDTesterDeviceNet](#) ()**

Initializes a new instance of the [CExternDTesterDeviceNet](#) class.

10.22.2.2 ~CExternDTesterDeviceNet() virtual [~CExternDTesterDeviceNet](#) () [virtual]**10.22.2.3 "!CExternDTesterDeviceNet() [!CExternDTesterDeviceNet](#) ()****10.22.3 Member Function Documentation****10.22.3.1 Read() array<uint8_t> ^ Read (int configString_Length)**

Reads the config string from the device.

Parameters

<i>configString_Length</i>	The maximal length of configString.
----------------------------	-------------------------------------

Returns

The config string.

10.22.3.2 Read2() `String ^ Read2 ()`

Reads the config string from the device.

Returns

The config string.

10.22.3.3 Write() `void Write (array< uint8_t >^ configString)`

Reads the config string from the device.

Parameters

<i>configString</i>	The config string.
---------------------	--------------------

10.22.3.4 Write2() `void Write2 (String^ configString)`

Reads the config string from the device.

Parameters

<i>configString</i>	The config string.
---------------------	--------------------

10.23 CFilterCoefficientsNet Class Reference**Public Member Functions**

- void [SetAFormat](#) (int vk, int nk, int pos, FilterCalculationDirectionEnumNet dir)
- void [SetBFormat](#) (int vk, int nk, int pos, FilterCalculationDirectionEnumNet dir)

- [CFilterCoefficientsNet](#) ()
- [CFilterCoefficientsNet](#) (double b0, double b1, double b2, double a1, double a2)
- [CFilterCoefficientsNet](#) (double b0, double b1, double a1)
- [CFilterCoefficientsNet](#) (array< double >^ b, array< double >^ a)
- [CFilterCoefficientsNet](#) (uint64_t b0, uint64_t b1, uint64_t b2, uint64_t a1, uint64_t a2)
- [~CFilterCoefficientsNet](#) ()
- virtual bool [Equals](#) ([CFilterCoefficientsNet](#)^ coefficients)

Properties

- array< double >^ [A](#) [get]
- array< double >^ [B](#) [get]
- double [B0](#) [get]
- double [B1](#) [get]
- double [B2](#) [get]
- double [A1](#) [get]
- double [A2](#) [get]
- uint64_t [UIntB0](#) [get]
- uint64_t [UIntB1](#) [get]
- uint64_t [UIntB2](#) [get]
- uint64_t [UIntA1](#) [get]
- uint64_t [UIntA2](#) [get]

10.23.1 Constructor & Destructor Documentation

10.23.1.1 CFilterCoefficientsNet() [1/5] [CFilterCoefficientsNet](#) ()

10.23.1.2 CFilterCoefficientsNet() [2/5] [CFilterCoefficientsNet](#) (
double b0,
double b1,
double b2,
double a1,
double a2)

10.23.1.3 CFilterCoefficientsNet() [3/5] [CFilterCoefficientsNet](#) (
double b0,
double b1,
double a1)

10.23.1.4 CFilterCoefficientsNet() [4/5] `CFilterCoefficientsNet` (
 array< double >^ *b*,
 array< double >^ *a*)

10.23.1.5 CFilterCoefficientsNet() [5/5] `CFilterCoefficientsNet` (
 uint64_t *b0*,
 uint64_t *b1*,
 uint64_t *b2*,
 uint64_t *a1*,
 uint64_t *a2*)

10.23.1.6 ~CFilterCoefficientsNet() `~CFilterCoefficientsNet` ()

10.23.2 Member Function Documentation

10.23.2.1 Equals() virtual bool Equals (
 `CFilterCoefficientsNet`^ *coefficients*) [virtual]

10.23.2.2 SetAFormat() void SetAFormat (
 int *vk*,
 int *nk*,
 int *pos*,
 FilterCalculationDirectionEnumNet *dir*)

10.23.2.3 SetBFormat() void SetBFormat (
 int *vk*,
 int *nk*,
 int *pos*,
 FilterCalculationDirectionEnumNet *dir*)

10.23.3 Property Documentation

10.23.3.1 A array< double>^ *A* [get]

10.23.3.2 A1 double A1 [get]

10.23.3.3 A2 double A2 [get]

10.23.3.4 B array< double>^ B [get]

10.23.3.5 B0 double B0 [get]

10.23.3.6 B1 double B1 [get]

10.23.3.7 B2 double B2 [get]

10.23.3.8 UIntA1 uint64_t UIntA1 [get]

10.23.3.9 UIntA2 uint64_t UIntA2 [get]


10.23.3.10 UIntB0 uint64_t UIntB0 [get]

10.23.3.11 UIntB1 uint64_t UIntB1 [get]

10.23.3.12 UIntB2 uint64_t UIntB2 [get]

10.24 CFilterConfigurationNet Class Reference

Inheritance diagram for CFilterConfigurationNet:

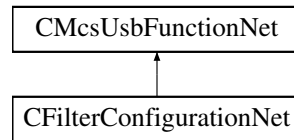

Public Member Functions

- [CFilterConfigurationNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetFilterParameter](#) (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, [CFilterCoefficientsNet](#)[^] Coefficients, [CFilterPropertyNet](#)[^] FilterProp)
- void [SetFilterParameter](#) (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, [CFilterCoefficientsNet](#)[^] CoefficientsSet1, [CFilterCoefficientsNet](#)[^] CoefficientsSet2, [CFilterPropertyNet](#)[^] FilterProp)
- void [SetFilterParameterPermanent](#) (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)
- void [EraseFilterParameterPermanent](#) (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)
- void [SetHighpassFilterEnable](#) (bool enable)
- bool [GetHighpassFilterEnable](#) ()
- void [ResetHighpassFilter](#) ()

Additional Inherited Members

10.24.1 Constructor & Destructor Documentation

10.24.1.1 CFilterConfigurationNet() [CFilterConfigurationNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.24.2 Member Function Documentation

10.24.2.1 EraseFilterParameterPermanent() void EraseFilterParameterPermanent (
DacqGroupChannelEnumNet GroupID,
uint32_t FilterNumber)

10.24.2.2 GetHighpassFilterEnable() bool GetHighpassFilterEnable ()

10.24.2.3 ResetHighpassFilter() void ResetHighpassFilter ()

10.24.2.4 SetFilterParameter() [1/2] void SetFilterParameter (
DacqGroupChannelEnumNet GroupID,
uint32_t FilterNumber,
CFilterCoefficientsNet^ Coefficients,
CFilterPropertyNet^ FilterProp)


10.24.2.5 SetFilterParameter() [2/2] void SetFilterParameter (
DacqGroupChannelEnumNet GroupID,
uint32_t FilterNumber,
CFilterCoefficientsNet^ CoefficientsSet1,
CFilterCoefficientsNet^ CoefficientsSet2,
CFilterPropertyNet^ FilterProp)

10.24.2.6 SetFilterParameterPermanent() void SetFilterParameterPermanent (
DacqGroupChannelEnumNet GroupID,
uint32_t FilterNumber)

10.24.2.7 SetHighpassFilterEnable() void SetHighpassFilterEnable (
bool enable)

10.25 CFilterConfigurationRegisterNet Class Reference

Inheritance diagram for CFilterConfigurationRegisterNet:

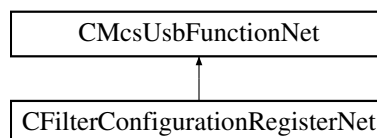

Public Member Functions

- CFilterConfigurationRegisterNet (CMcsUsbNet^ mcsusb)
- void SetFilterParameter (uint32_t FilterCoefRegBase, CFilterCoefficientsNet^ Coefficients, uint32_t FilterInfoRegBase, CFilterPropertyNet^ FilterProp)
- void SetFilterParameter (uint32_t FilterCoefSet1RegBase, CFilterCoefficientsNet^ CoefficientsSet1, uint32_t FilterCoefSet2RegBase, CFilterCoefficientsNet^ CoefficientsSet2, uint32_t FilterInfoRegBase, CFilterPropertyNet^ FilterProp)
- void SetFilterParameterPermanent (uint32_t FilterCoefRegBase, uint32_t FilterCoefDmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void SetFilterParameterPermanent (uint32_t FilterCoefSet1RegBase, uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2RegBase, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefDmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)

Additional Inherited Members

10.25.1 Constructor & Destructor Documentation

10.25.1.1 CFilterConfigurationRegisterNet() [CFilterConfigurationRegisterNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

10.25.2 Member Function Documentation

10.25.2.1 EraseFilterParameterPermanent() [1/2] void EraseFilterParameterPermanent (
 uint32_t *FilterCoefDmaReg*,
 uint32_t *FilterInfoDmaReg*,
 uint32_t *EEPROMBase*,
 uint32_t *EEPROMSize*)

10.25.2.2 EraseFilterParameterPermanent() [2/2] void EraseFilterParameterPermanent (
 uint32_t *FilterCoefSet1DmaReg*,
 uint32_t *FilterCoefSet2DmaReg*,
 uint32_t *FilterInfoDmaReg*,
 uint32_t *EEPROMBase*,
 uint32_t *EEPROMSize*)

10.25.2.3 SetFilterParameter() [1/2] void SetFilterParameter (
 uint32_t *FilterCoefRegBase*,
[CFilterCoefficientsNet](#)[^] *Coefficients*,
 uint32_t *FilterInfoRegBase*,
[CFilterPropertyNet](#)[^] *FilterProp*)

10.25.2.4 SetFilterParameter() [2/2] void SetFilterParameter (
 uint32_t *FilterCoefSet1RegBase*,
[CFilterCoefficientsNet](#)[^] *CoefficientsSet1*,
 uint32_t *FilterCoefSet2RegBase*,
[CFilterCoefficientsNet](#)[^] *CoefficientsSet2*,
 uint32_t *FilterInfoRegBase*,
[CFilterPropertyNet](#)[^] *FilterProp*)

10.25.2.5 SetFilterParameterPermanent() [1/2] `void SetFilterParameterPermanent (`
`uint32_t FilterCoefRegBase,`
`uint32_t FilterCoefDmaReg,`
`uint32_t FilterInfoRegBase,`
`uint32_t FilterInfoDmaReg,`
`uint32_t EEPROMBase,`
`uint32_t EEPROMSize)`

10.25.2.6 SetFilterParameterPermanent() [2/2] `void SetFilterParameterPermanent (`
`uint32_t FilterCoefSet1RegBase,`
`uint32_t FilterCoefSet1DmaReg,`
`uint32_t FilterCoefSet2RegBase,`
`uint32_t FilterCoefSet2DmaReg,`
`uint32_t FilterInfoRegBase,`
`uint32_t FilterInfoDmaReg,`
`uint32_t EEPROMBase,`
`uint32_t EEPROMSize)`

10.26 CFilterPropertyNet Class Reference

Public Member Functions

- [CFilterPropertyNet](#) (uint32_t CornerFrequenzymHz, uint32_t [Order](#), FilterBandEnumNet [FilterBand](#), FilterFamilyEnumNet [FilterFamily](#), FilterTypeEnumNet [FilterType](#), bool Active)
- [~CFilterPropertyNet](#) ()
- virtual System::String ^ [ToString](#) () override

Properties

- uint32_t [CornerFrequencymHz](#) [get]
- double [CornerFrequency](#) [get]
- uint32_t [Order](#) [get]
- FilterBandEnumNet [FilterBand](#) [get]
- FilterFamilyEnumNet [FilterFamily](#) [get]
- FilterTypeEnumNet [FilterType](#) [get]
- bool [FilterActive](#) [get]

10.26.1 Constructor & Destructor Documentation

10.26.1.1 CFilterPropertyNet() `CFilterPropertyNet (`
`uint32_t CornerFrequenzymHz,`
`uint32_t Order,`
`FilterBandEnumNet FilterBand,`
`FilterFamilyEnumNet FilterFamily,`
`FilterTypeEnumNet FilterType,`
`bool Active)`

10.26.1.2 `~CFilterPropertyNet()` `~CFilterPropertyNet ()`

10.26.2 Member Function Documentation

10.26.2.1 `ToString()` `virtual System::String ^ ToString () [override], [virtual]`

10.26.3 Property Documentation

10.26.3.1 `CornerFrequency` `double CornerFrequency [get]`

10.26.3.2 `CornerFrequencymHz` `uint32_t CornerFrequencymHz [get]`

10.26.3.3 `FilterActive` `bool FilterActive [get]`

10.26.3.4 `FilterBand` `FilterBandEnumNet FilterBand [get]`

10.26.3.5 `FilterFamily` `FilterFamilyEnumNet FilterFamily [get]`


10.26.3.6 `FilterType` `FilterTypeEnumNet FilterType [get]`

10.26.3.7 `Order` `uint32_t Order [get]`

10.27 CFluidControlDeviceNet Class Reference

[CFluidControlDeviceNet](#) is the class to control MCS FluidControl (FCB and FCX) device.

Inheritance diagram for CFluidControlDeviceNet:

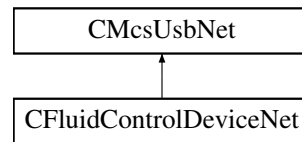

Public Member Functions

- [CFluidControlDeviceNet](#) ()
Initialize a new instance of the [CFluidControlDeviceNet](#) class.
- [~CFluidControlDeviceNet](#) ()
Default destructor.
- void [SetValve](#) (unsigned int value)
Open or Close valves.
- void [SetSingleValve](#) (unsigned short valve, unsigned short onoff)
Opens or Closes a valve.
- void [SetDigout](#) (unsigned int value)
Define the pattern on the Digital Output.
- void [SetPWM](#) (unsigned int channel, unsigned int value)
Sets the duty cycle of the PWM output.
- void [CalibrateThermocouple](#) (unsigned int channel)
Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.
- void [SetThermocoupleNanovoltPerKelvin](#) (unsigned int channel, unsigned int value)
Sets the proportional constant for the Thermocouple.
- unsigned int [GetValve](#) ()
Gets the state of the valves.
- unsigned short [GetSingleValve](#) (unsigned short valve)
Gets the state of a valve.
- unsigned int [GetDigout](#) ()
Gets the state of the digital output.
- unsigned int [GetPWM](#) (unsigned int channel)
Gets the state of the PWM output.
- unsigned int [GetAdc](#) (unsigned int channel)
Reads an ADC Value.
- unsigned int [GetDigin](#) ()
Reads the digital input.
- int [GetThermocoupleTemperature](#) (unsigned int channel)
Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.
- int [GetReferenceTemperature](#) (unsigned int channel)
Reads the reference temperature for the Thermocouple.
- unsigned int [GetThermocoupleCalibration](#) (unsigned int channel)
Gets the calibration constant for the Thermocouple ADC.
- unsigned int [GetThermocoupleNanovoltPerKelvin](#) (unsigned int channel)
Reads the proportional constant for the Thermocouple.

Properties

- [CMcsBus_VoltageModeNet](#)[^] [McsBus_VoltageMode](#) [get]

Additional Inherited Members

10.27.1 Detailed Description

[CFluidControlDeviceNet](#) is the class to control MCS FluidControl (FCB and FCX) device.

10.27.2 Constructor & Destructor Documentation

10.27.2.1 CFluidControlDeviceNet() [CFluidControlDeviceNet](#) ()

Initialize a new instance of the [CFluidControlDeviceNet](#) class.

10.27.2.2 ~CFluidControlDeviceNet() [~CFluidControlDeviceNet](#) ()

Default destructor.

10.27.3 Member Function Documentation

10.27.3.1 CalibrateThermocouple() [void CalibrateThermocouple](#) ([unsigned int channel](#))

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

10.27.3.2 GetAdc() [unsigned int GetAdc](#) ([unsigned int channel](#))

Reads an ADC Value.

Parameters

<i>channel</i>	The ADC channel number to query.
----------------	----------------------------------

Returns

The current ADC value.

10.27.3.3 GetDigin() `unsigned int GetDigin ()`

Reads the digital input.

Returns

The bit pattern of the state of the digital inputs.

10.27.3.4 GetDigout() `unsigned int GetDigout ()`

Gets the state of the digital output.

Returns

The current state of the digital outputs as a bit pattern.

10.27.3.5 GetPWM() `unsigned int GetPWM (`
`unsigned int channel)`

Gets the state of the PWM output.

Returns

The current state of the PWM outputs duty cycle in permille.

10.27.3.6 GetReferenceTemperature() `int GetReferenceTemperature (`
`unsigned int channel)`

Reads the reference temperature for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The temperature from the Thermocouple in 1/100 °C.

10.27.3.7 GetSingleValve() unsigned short GetSingleValve (
unsigned short valve)

Gets the state of a valve.

Parameters

<i>valve</i>	number of valve
--------------	-----------------

Returns

state of the valve

10.27.3.8 GetThermocoupleCalibration() unsigned int GetThermocoupleCalibration (
unsigned int channel)

Gets the calibration constant for the Thermocouple ADC.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The calibration constant for the Thermocouple ADC.

10.27.3.9 GetThermocoupleNanovoltPerKelvin() unsigned int GetThermocoupleNanovoltPerKelvin (
unsigned int channel)

Reads the proportional constant for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

10.27.3.10 GetThermocoupleTemperature() `int GetThermocoupleTemperature (unsigned int channel)`

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The temperature difference between both Thermocouple junctions in 1/100 °C.

10.27.3.11 GetValve() `unsigned int GetValve ()`

Gets the state of the valves.

Returns

The current state of the valves as a bit pattern.

10.27.3.12 SetDigout() `void SetDigout (unsigned int value)`

Define the pattern on the Digital Output.

Parameters

<i>value</i>	bit pattern on the digital output.
--------------	------------------------------------

10.27.3.13 SetPWM() `void SetPWM (unsigned int channel, unsigned int value)`

Sets the duty cycle of the PWM output.

Parameters

<i>channel</i>	PWM channel number.
<i>value</i>	duty cycle of the PWM output in permille.

10.27.3.14 SetSingleValve() `void SetSingleValve (`
 `unsigned short valve,`
 `unsigned short onoff)`

Opens or Closes a valve.

Parameters

<i>valve</i>	number of valve to be changed.
--------------	--------------------------------

Parameters

<i>onoff</i>	open or close the valve.
--------------	--------------------------

10.27.3.15 SetThermocoupleNanovoltPerKelvin() `void SetThermocoupleNanovoltPerKelvin (`
 `unsigned int channel,`
 `unsigned int value)`

Sets the proportinal constant for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
<i>value</i>	proportinal constant for the Thermocouple in Nanovolt per Kelvin.

10.27.3.16 SetValve() `void SetValve (`
 `unsigned int value)`

Open or Close valves.

Parameters

<i>value</i>	bit pattern of valves which should be open.
--------------	---


10.27.4 Property Documentation

10.27.4.1 McsBus_VoltageMode [CMcsBus_VoltageModeNet](#)[^] McsBus_VoltageMode [get]

10.28 CFYIDeviceNet Class Reference

[CFYIDeviceNet](#) is the class to control the MCS FYI device

Inheritance diagram for CFYIDeviceNet:

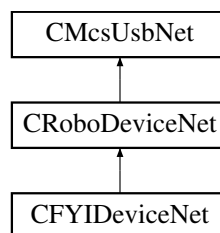

Public Member Functions

- [CFYIDeviceNet](#) (void)

Properties

- [CRobo_FYITemp_FunctionNet](#)[^] FYITemp [get]
- [CRobo_FYIProgram_FunctionNet](#)[^] FYIProgram [get]
- [CMcsBus_SensorNet](#)[^] Sensor [get]

Additional Inherited Members

10.28.1 Detailed Description

[CFYIDeviceNet](#) is the class to control the MCS FYI device

10.28.2 Constructor & Destructor Documentation

10.28.2.1 CFYIDeviceNet() [CFYIDeviceNet](#) (
void)

10.28.3 Property Documentation

10.28.3.1 FYIProgram `CRobo_FYIProgram_FunctionNet^ FYIProgram [get]`


10.28.3.2 FYITemp `CRobo_FYITemp_FunctionNet^ FYITemp [get]`

10.28.3.3 Sensor `CMcsBus_SensorNet^ Sensor [get]`

10.29 CGenericDevelopDeviceNet Class Reference

`CGenericDevelopDeviceNet` is the class to use during development of a new device.

Inheritance diagram for `CGenericDevelopDeviceNet`:

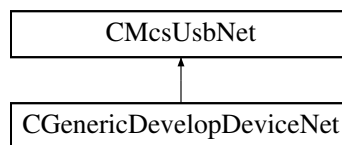

Public Member Functions

- `CGenericDevelopDeviceNet` (void)
Initialize a new instance of the `CGenericDevelopDeviceNet` class.
- `~CGenericDevelopDeviceNet` (void)
- void `SetValue` (uint16_t value, uint16_t index)
Sets .

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

- `template<typename C >`
void `SetBuffer` (uint16_t value, uint16_t index, array< C >^ buffer)

- void [SetUByteBuffer](#) (uint16_t value, uint16_t index, array< unsigned char >^ buffer)

Sends an array of type unsigned char to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetByteBuffer](#) (uint16_t value, uint16_t Index, array< char >^ buffer)

Sends an array of type char to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

Index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetUShortBuffer](#) (uint16_t value, uint16_t index, array< unsigned short >^ buffer)

Sends an array of type unsigned short to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetShortBuffer](#) (uint16_t value, uint16_t index, array< short >^ buffer)
Sends an array of type short to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetUIntBuffer](#) (uint16_t value, uint16_t index, array< unsigned int >^ buffer)
Sends an array of unsigned int to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetIntBuffer](#) (uint16_t value, uint16_t index, array< int >^ buffer)
Sends an array of type int to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- `template<typename C >`
`array< C > ^ GetBuffer (uint16_t value, uint16_t index, int size)`
- `array< unsigned char > ^ GetUByteBuffer (uint16_t value, uint16_t index, int size)`
Gets an array of type unsigned char from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- `array< char > ^ GetByteBuffer (uint16_t value, uint16_t index, int size)`
Gets an array of type char from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< unsigned short > ^ [GetUShortBuffer](#) (uint16_t value, uint16_t index, int size)
Gets an array of type unsigned short from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< short > ^ [GetShortBuffer](#) (uint16_t value, uint16_t index, int size)
Gets an array of type short from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< unsigned int > ^ [GetUIntBuffer](#) (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned int from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< int > ^ [GetIntBuffer](#) (uint16_t value, uint16_t index, int size)

Gets an array of type int from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

*Returns**The array of data from the device.*

- IntPtr [OpenPipe](#) (uint8_t endpointAddress)

*Open a Pipe to an USB Endpoint.**Parameters*

endpointAddress	<i>The Endpoint Number.</i>
-----------------	-----------------------------

*Returns**A handle to the endpoint.*

- void [ClosePipe](#) (IntPtr pipeHandle)

*Close a Pipe to an USB Endpoint.**Parameters*

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

- void [ResetPipe](#) (IntPtr pipeHandle)

*Reset a Pipe to an USB Endpoint.**Parameters*

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

- template<typename C >
array< C > ^ [ReadPipe](#) (IntPtr pipeHandle, uint32_t size)

*Read data from an USB Endpoint.**Parameters*

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

Parameters

size	<i>Number of elements to read.</i>
------	------------------------------------

*Returns**An array of data read.*

- template<typename C >
void [WritePipe](#) (IntPtr pipeHandle, array< C >^ buffer)

*Write data to an USB Endpoint.**Parameters*

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

Parameters

<code>buffer</code>	<i>An array of data to write.</i>
---------------------	-----------------------------------

Additional Inherited Members**10.29.1 Detailed Description**

[CGenericDevelopDeviceNet](#) is the class to use during development of a new device.

10.29.2 Constructor & Destructor Documentation

10.29.2.1 CGenericDevelopDeviceNet() [CGenericDevelopDeviceNet](#) (
 void)

Initialize a new instance of the [CGenericDevelopDeviceNet](#) class.

10.29.2.2 ~CGenericDevelopDeviceNet() [~CGenericDevelopDeviceNet](#) (
 void)

10.29.3 Member Function Documentation

10.29.3.1 ClosePipe() void ClosePipe (
 IntPtr *pipeHandle*)

Close a Pipe to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

10.29.3.2 GetBuffer() `array<C> ^ GetBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

10.29.3.3 GetByteBuffer() `array<char> ^ GetByteBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type char from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.4 GetIntBuffer() `array<int> ^ GetIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type int from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.5 GetShortBuffer() `array<short> ^ GetShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type short from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.6 GetUByteBuffer() `array<unsigned char> ^ GetUByteBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type unsigned char from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.7 GetUIntBuffer() `array<unsigned int> ^ GetUIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type unsigned int from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.8 GetUShortBuffer() `array<unsigned short> ^ GetUShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type unsigned short from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

10.29.3.9 OpenPipe() `IntPtr OpenPipe (`
`uint8_t endpointAddress)`

Open a Pipe to an USB Endpoint.

Parameters

<i>endpointAddress</i>	The Endpoint Number.
------------------------	----------------------

Returns

A handle to the endpoint.

10.29.3.10 ReadPipe() `array<C> ^ ReadPipe (`
`IntPtr pipeHandle,`
`uint32_t size)`

Read data from an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

Parameters

<i>size</i>	Number of elements to read.
-------------	-----------------------------

Returns

An array of data read.

10.29.3.11 ResetPipe() `void ResetPipe (`
 `IntPtr pipeHandle)`

Reset a Pipe to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

10.29.3.12 SetBuffer() `void SetBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< C >^ buffer)`

10.29.3.13 SetByteBuffer() `void SetByteBuffer (`
 `uint16_t value,`
 `uint16_t Index,`
 `array< char >^ buffer)`

Sends an array of type char to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>Index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.14 SetIntBuffer() `void SetIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< int >^ buffer)`

Sends an array of type int to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.15 SetShortBuffer() `void SetShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< short >^ buffer)`

Sends an array of type short to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.16 SetUByteBuffer() void SetUByteBuffer (
 uint16_t *value*,
 uint16_t *index*,
 array< unsigned char >^ *buffer*)

Sends an array of type unsigned char to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.17 SetUIntBuffer() void SetUIntBuffer (
 uint16_t *value*,
 uint16_t *index*,
 array< unsigned int >^ *buffer*)

Sends an array of unsigned int to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.18 SetUShortBuffer() `void SetUShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< unsigned short >^ buffer)`

Sends an array of type unsigned short to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

10.29.3.19 SetValue() `void SetValue (`
 `uint16_t value,`
 `uint16_t index)`

Sets .

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

10.29.3.20 WritePipe() `void WritePipe (`
 `IntPtr pipeHandle,`
 `array< C >^ buffer)`

Write data to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------


Parameters

<i>buffer</i>	An array of data to write.
---------------	----------------------------

10.30 CGilsonDeviceNet Class Reference

[CGilsonDeviceNet](#) is the class to control a Gilson device.

Inheritance diagram for CGilsonDeviceNet:

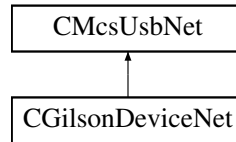

Public Member Functions

- [CGilsonDeviceNet](#) (void)
Initialize a new instance of the [CGilsonDeviceNet](#) class.
- [~CGilsonDeviceNet](#) (void)
- void [ConnectSlave](#) (byte ID)
- void [SendImmediate](#) (wchar_t command)
- String ^ [SendImmediateGetResponse](#) (wchar_t command)
- void [SendBuffered](#) (String^ command)
- String ^ [GetLastAnswer](#) ()

Protected Attributes

- CGilsonDevice * [m_pGilsonDevice](#)

Additional Inherited Members

10.30.1 Detailed Description

[CGilsonDeviceNet](#) is the class to control a Gilson device.

10.30.2 Constructor & Destructor Documentation

10.30.2.1 [CGilsonDeviceNet\(\)](#) [CGilsonDeviceNet](#) (
void)

Initialize a new instance of the [CGilsonDeviceNet](#) class.

10.30.2.2 [~CGilsonDeviceNet\(\)](#) [~CGilsonDeviceNet](#) (
void)

10.30.3 Member Function Documentation

10.30.3.1 ConnectSlave() `void ConnectSlave (`
`byte ID)`

10.30.3.2 GetLastAnswer() `String ^ GetLastAnswer ()`

10.30.3.3 SendBuffered() `void SendBuffered (`
`String^ command)`

10.30.3.4 SendImmediate() `void SendImmediate (`
`wchar_t command)`

10.30.3.5 SendImmediateGetResponse() `String ^ SendImmediateGetResponse (`
`wchar_t command)`


10.30.4 Member Data Documentation

10.30.4.1 m_pGilsonDevice `CGilsonDevice* m_pGilsonDevice [protected]`

10.31 CHiClampDeviceNet Class Reference

[CHiClampDeviceNet](#) is the to control the MCS HiClamp device

Inheritance diagram for CHiClampDeviceNet:

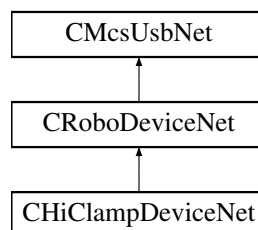

Public Member Functions

- [CHiClampDeviceNet](#) (void)

Properties

- [CRoboDacqNet](#)[^] [RoboDacq](#) [get]

Additional Inherited Members

10.31.1 Detailed Description

[CHiClampDeviceNet](#) is the to control the MCS HiClamp device

10.31.2 Constructor & Destructor Documentation


10.31.2.1 CHiClampDeviceNet() [CHiClampDeviceNet](#) (
void)

10.31.3 Property Documentation

10.31.3.1 RoboDacq [CRoboDacqNet](#)[^] [RoboDacq](#) [get]

10.32 CHLADacqNet Class Reference

Inheritance diagram for CHLADacqNet:

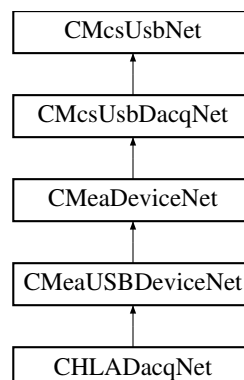

Public Member Functions

- [CHLADacqNet](#) (void)

Additional Inherited Members


10.32.1 Constructor & Destructor Documentation

10.32.1.1 CHLADacqNet() `CHLADacqNet (void)`

10.33 CHLADeviceNet Class Reference

[CHLADeviceNet](#) is the to control the MCS HLA device

Inheritance diagram for CHLADeviceNet:

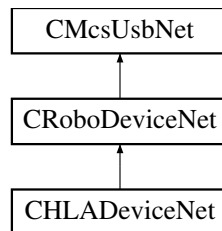

Public Member Functions

- [CHLADeviceNet](#) (void)

Properties

- [CHLADacqNet](#)[^] [HLADacq](#) [get]
- [CSerialPortNet](#)[^] [SerialPort](#) [get]

Additional Inherited Members

10.33.1 Detailed Description

[CHLADeviceNet](#) is the to control the MCS HLA device

10.33.2 Constructor & Destructor Documentation

10.33.2.1 CHLADeviceNet() `CHLADeviceNet (void)`

10.33.3 Property Documentation

10.33.3.1 HLADataq [CHLADacqNet](#)[^] HLADataq [get]

10.33.3.2 SerialPort [CSerialPortNet](#)[^] SerialPort [get]

10.34 CMcsUsbDacqNet::CHWInfo Class Reference

Class to provide hardware information about the device.

Classes

- class [CVoltageRangeInfoNet](#)

Public Member Functions

- [CHWInfo](#) ([CMcsUsbDacqNet](#)[^] device)
- virtual uint32_t [GetNumberOfHWADCCChannels](#) ([System::Runtime::InteropServices::Out]int% numberOfChannels)
Get the number of analog channels the device supports.
- virtual uint32_t [GetNumberOfHWDigitalChannels](#) ([System::Runtime::InteropServices::Out]int% numberOfChannels)
Get the number of digital channels the device supports.
- virtual bool [IsDigitalChannelDedicated](#) ()
Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.
- virtual uint32_t [GetAvailableSampleRates](#) ([System::Runtime::InteropServices::Out]System::Collections::Generic::List< int32_t > ^% sampleRates)
- virtual uint32_t [GetAvailableVoltageRangesInMicroVolt](#) ([System::Runtime::InteropServices::Out]System::Collections::Generic::List< int32_t > ^% voltageRanges)
- virtual uint32_t [GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt](#) ([System::Runtime::InteropServices::Out]System::Collections::Generic::List< [CVoltageRangeInfoNet](#)[^] > ^% voltageRanges)

10.34.1 Detailed Description

Class to provide hardware information about the device.

10.34.2 Constructor & Destructor Documentation

10.34.2.1 CHWInfo() [CHWInfo](#) (
[CMcsUsbDacqNet](#)[^] *device*)

10.34.3 Member Function Documentation

10.34.3.1 GetAvailableSampleRates() [virtual uint32_t GetAvailableSampleRates](#) (
[System::Runtime::InteropServices::Out] [System::Collections::Generic::List](#)< [int32](#)↔
_t >^% *sampleRates*) [virtual]

10.34.3.2 GetAvailableVoltageRangesInMicroVolt() [virtual uint32_t GetAvailableVoltageRangesIn](#)↔
[MicroVolt](#) (
[System::Runtime::InteropServices::Out] [System::Collections::Generic::List](#)< [int32](#)↔
_t >^% *voltageRanges*) [virtual]

10.34.3.3 GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt() [virtual uint32_t GetAvailable](#)↔
[VoltageRangesInMicroVoltAndStringsInMilliVolt](#) (
[System::Runtime::InteropServices::Out] [System::Collections::Generic::List](#)< [CVoltageRangeInfoNet](#)
voltageRanges) [virtual]

10.34.3.4 GetNumberOfHWADCCChannels() [virtual uint32_t GetNumberOfHWADCCChannels](#) (
[System::Runtime::InteropServices::Out] [int](#)% *numberOfChannels*) [virtual]

Get the number of analog channels the device supports.

Parameters

<i>numberOfChannels</i>	Number of analog channels the device supports.
-------------------------	--

Returns

Error Status. 0 on success.

10.34.3.5 GetNumberOfHWDigitalChannels() [virtual uint32_t GetNumberOfHWDigitalChannels](#) (
[System::Runtime::InteropServices::Out] [int](#)% *numberOfChannels*) [virtual]

Get the number of digital channels the device supports.

Parameters

<i>numberOfChannels</i>	Number of digital channels the device supports.
-------------------------	---

Returns

Error Status. 0 on success.

10.34.3.6 IsDigitalChannelDedicated() `virtual bool IsDigitalChannelDedicated () [virtual]`


Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

Returns

false when the digital channel replaces an analog channel when enabled, true when the digital channels is appended to the analog channels when enabled.

10.35 CIntanMea_FunctionNet Class Reference

Inheritance diagram for CIntanMea_FunctionNet:

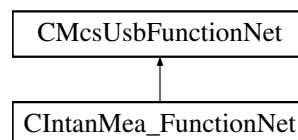

Public Member Functions

- [CIntanMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] intalMea_FunctionPointerContainer)
- [CIntanMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- [int](#) [GetUpperFrequencyByIndex](#) (unsigned short index)
- [int](#) [GetLowerFrequencyByIndex](#) (unsigned short index)
- [int64_t](#) [GetDSPHighPassByIndex](#) (unsigned short index)
- [int](#) [GetIntanRegister](#) (unsigned short chip, unsigned short registernumber)
- [int](#) [GetImpedanceResult](#) (unsigned short channel)
- [void](#) [SetBandwidthByIndex](#) (int upper_index, int lower_index)
- [void](#) [SetDSPHighPassByIndex](#) (int index)
- [void](#) [AmplifierSettle](#) ()
- [void](#) [SetIntanRegister](#) (unsigned short register_number, int value)
- [void](#) [SetDiagnosticMode](#) (unsigned char onoff)
- [void](#) [BeginImpedanceCheck](#) (array< int >[^] config_values)

Additional Inherited Members

10.35.1 Constructor & Destructor Documentation

10.35.1.1 CIntanMea_FunctionNet() [1/2] `CIntanMea_FunctionNet (`
`CMcsUsbNet^ mcsusb,`
`CMcsUsbFunctionPointerContainer^ intalMea_FunctionPointerContainer)`

10.35.1.2 CIntanMea_FunctionNet() [2/2] `CIntanMea_FunctionNet (`
`CMcsUsbNet^ mcsusb)`

10.35.2 Member Function Documentation

10.35.2.1 AmplifierSettle() `void AmplifierSettle ()`

10.35.2.2 BeginImpedanceCheck() `void BeginImpedanceCheck (`
`array< int >^ config_values)`

10.35.2.3 GetDSPHighPassByIndex() `int64_t GetDSPHighPassByIndex (`
`unsigned short index)`

10.35.2.4 GetImpedanceResult() `int GetImpedanceResult (`
`unsigned short channel)`

10.35.2.5 GetIntanRegister() `int GetIntanRegister (`
`unsigned short chip,`
`unsigned short registernumber)`

10.35.2.6 GetLowerFrequencyByIndex() `int GetLowerFrequencyByIndex (`
`unsigned short index)`

10.35.2.7 GetUpperFrequencyByIndex() `int GetUpperFrequencyByIndex (`
`unsigned short index)`

10.35.2.8 SetBandwidthByIndex() `void SetBandwidthByIndex (`
`int upper_index,`
`int lower_index)`

10.35.2.9 SetDiagnosticMode() `void SetDiagnosticMode (`
`unsigned char onoff)`


10.35.2.10 SetDSPHighPassByIndex() `void SetDSPHighPassByIndex (`
`int index)`

10.35.2.11 SetIntanRegister() `void SetIntanRegister (`
`unsigned short register_number,`
`int value)`

10.36 CInterfaceboardFunctionNet Class Reference

[CInterfaceboardFunctionNet](#) is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboardFunctionNet:

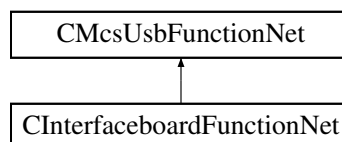

Public Member Functions

- [CInterfaceboardFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pInterfaceboardFunctionPointerContainer)
Initializes a new instance of the [CInterfaceboardFunctionNet](#) class.
- [CInterfaceboardFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual [~CInterfaceboardFunctionNet](#) ()
- [!CInterfaceboardFunctionNet](#) ()
- void [SetCardinalDacqSamplerate](#) (uint32_t samplerate)
Sets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz
- uint32_t [GetCardinalDacqSamplerate](#) ()
Gets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz
- void [SetCardinalStgOutputrate](#) (uint32_t outputrate)
Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz
- uint32_t [GetCardinalStgOutputrate](#) ()
Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Additional Inherited Members

10.36.1 Detailed Description

[CInterfaceboardFunctionNet](#) is the class to control the Interfaceboard

10.36.2 Constructor & Destructor Documentation

10.36.2.1 CInterfaceboardFunctionNet() [1/2] [CInterfaceboardFunctionNet](#) (
[CMcsUsbNet](#)^ mcsusb,
[CMcsUsbFunctionPointerContainer](#)^ pInterfaceboardFunctionPointerContainer)

Initializes a new instance of the [CInterfaceboardFunctionNet](#) class.

10.36.2.2 CInterfaceboardFunctionNet() [2/2] [CInterfaceboardFunctionNet](#) (
[CMcsUsbNet](#)^ mcsusb)

10.36.2.3 ~CInterfaceboardFunctionNet() virtual [~CInterfaceboardFunctionNet](#) () [virtual]

10.36.2.4 !CInterfaceboardFunctionNet() [!CInterfaceboardFunctionNet](#) ()

10.36.3 Member Function Documentation

10.36.3.1 GetCardinalDacqSamplerate() `uint32_t GetCardinalDacqSamplerate ()`

Gets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz

Returns

The samplerate in Hz.

10.36.3.2 GetCardinalStgOutputrate() `uint32_t GetCardinalStgOutputrate ()`

Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Returns

The output rate in Hz.

10.36.3.3 SetCardinalDacqSamplerate() `void SetCardinalDacqSamplerate (uint32_t samplerate)`

Sets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz

Parameters

<i>samplerate</i>	The samplerate in Hz.
-------------------	-----------------------

10.36.3.4 SetCardinalStgOutputrate() `void SetCardinalStgOutputrate (uint32_t outputrate)`

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz


Parameters

<i>outputrate</i>	The output rate in Hz.
-------------------	------------------------

10.37 CLIH3DeviceNet Class Reference

[CLIH3DeviceNet](#) is the class to access the HEKA LIH3 device.

Inheritance diagram for CLIH3DeviceNet:

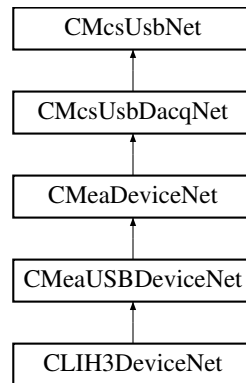

Public Member Functions

- [CLIH3DeviceNet](#) ()
Initializes a new instance of the [CLIH3DeviceNet](#) class.
- virtual [~CLIH3DeviceNet](#) ()
- [!CLIH3DeviceNet](#) ()
- void [DummyCommand](#) (uint32_t dummyParameter)
Dummy command to show how to use the DLL.
- void [SetEepromPage](#) (uint32_t EepromStartAddress, array< int8_t >^ EepromData, LIH30_EPC10_Bus_EnumNet epc10bus)
Writes into Eeprom on the EPC10 EEPROM
- array< int8_t >^ [GetEepromPage](#) (uint32_t EepromStartAddress, int EepromData_Length, LIH30_EPC10_Bus_EnumNet epc10bus)
Reads the requested amount of Eeprom byte from the EPC10 EEPROM
- void [SetSampleInterval](#) (uint32_t SampleInterval)
Sets the Sample Interval for the DACQ and Stimulation
- uint32_t [GetSampleInterval](#) ()
Gets the Sample Interval for the DACQ and Stimulation
- void [SetADCOffset](#) (LIH30_ADC_Channel_EnumNet AdcChannel, uint32_t Offset)
Sets the ADC offset of the DACQ for a single channel
- uint32_t [GetADCOffset](#) (LIH30_ADC_Channel_EnumNet AdcChannel)
Gets the ADC offset of the DACQ for a single channel
- uint32_t [ReadClipping](#) (LIH30_EPC10_Bus_EnumNet epc10bus)
Gets the clipping information
- void [SetDigOutState](#) (uint16_t DigOutState)
Writes to the LIH30 digital output
- uint16_t [GetDigInState](#) ()
Reads from the LIH30 digital input
- void [SendCommand](#) (LIH30_EPC10_Bus_EnumNet epc10bus, uint16_t Command)
Send command to the EPC10
- uint16_t [GetDacqRunStatus](#) ()
Gets the data acquisition running status
- void [SetDacUseldleValue](#) (uint32_t DacChannel, bool Useldle)
Sets if the DAC Idle value is used after stimulation
- bool [GetDacUseldleValue](#) (uint32_t DacChannel)
Gets if the DAC Idle value is used after stimulation

- void [SetDacIdleValue](#) (uint32_t DacChannel, int32_t IdleValue)
Sets the DAC Idle value
- int32_t [GetDacIdleValue](#) (uint32_t DacChannel)
Gets the DAC Idle value
- void [EnableUserTrigger](#) (bool enable)
Enables the User Trigger
- bool [IsUserTriggerEnabled](#) ()
Is the User Trigger enabled

Properties

- [CStimulusFunctionNet](#)[^] [StimulusFunction](#) [get]

Additional Inherited Members

10.37.1 Detailed Description

[CLIH3DeviceNet](#) is the class to access the HEKA LIH3 device.

10.37.2 Constructor & Destructor Documentation

10.37.2.1 CLIH3DeviceNet() [CLIH3DeviceNet](#) ()

Initializes a new instance of the [CLIH3DeviceNet](#) class.

10.37.2.2 ~CLIH3DeviceNet() [virtual](#) [~CLIH3DeviceNet](#) () [virtual]

10.37.2.3 "!CLIH3DeviceNet() [!CLIH3DeviceNet](#) ()

10.37.3 Member Function Documentation

10.37.3.1 DummyCommand() [void](#) [DummyCommand](#) (uint32_t *dummyParameter*)

Dummy command to show how to use the DLL.

Parameters

<i>dummyParameter</i>	parameter to send to the device
-----------------------	---------------------------------

10.37.3.2 EnableUserTrigger() `void EnableUserTrigger (`
 `bool enable)`

Enables the User Trigger

Parameters

<i>enable</i>	Enable
---------------	--------

10.37.3.3 GetADCOffset() `uint32_t GetADCOffset (`
 `LIH30_ADC_Channel_EnumNet AdcChannel)`

Gets the ADC offset of the DACQ for a single channel

Parameters

<i>AdcChannel</i>	The ADC channel
-------------------	-----------------

Returns

The offset for the given channel number

10.37.3.4 GetDacIdleValue() `int32_t GetDacIdleValue (`
 `uint32_t DacChannel)`

Gets the DAC Idle value

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

Returns

The idle value

10.37.3.5 GetDacqRunStatus() `uint16_t GetDacqRunStatus ()`

Gets the data acquisition running status

Returns

The status (1: running / 0: stopped)

10.37.3.6 GetDacUseIdleValue() `bool GetDacUseIdleValue (
uint32_t DacChannel)`

Gets if the DAC Idle value is used after stimulation

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

Returns

Use idle value

10.37.3.7 GetDigInState() `uint16_t GetDigInState ()`

Reads from the LIH30 digital input

Returns

The bit mask defining the digital input state

10.37.3.8 GetEepromPage() `array<int8_t> ^ GetEepromPage (
uint32_t EepromStartAddress,
int EepromData_Length,
LIH30_EPC10_Bus_EnumNet epc10bus)`

Reads the requested amount of Eeprom byte from the EPC10 EEPROM

Parameters

<i>EepromStartAddress</i>	start address of memory area to read from
<i>EepromData_Length</i>	The maximal length of EepromData.
<i>epc10bus</i>	The EPC10 bus

Returns

pointer to internal memory for the requested amount of data

10.37.3.9 GetSampleInterval() `uint32_t GetSampleInterval ()`

Gets the Sample Interval for the DACQ and Stimulation

Returns

Sample Interval configured on the device

10.37.3.10 IsUserTriggerEnabled() `bool IsUserTriggerEnabled ()`

Is the User Trigger enabled

Returns

Enabled

10.37.3.11 ReadClipping() `uint32_t ReadClipping (`
`LIH30_EPC10_Bus_EnumNet epc10bus)`

Gets the clipping information

Parameters

<i>epc10bus</i>	The EPC10 bus
-----------------	---------------

Returns

The clipping value

10.37.3.12 SendCommand() `void SendCommand (`
`LIH30_EPC10_Bus_EnumNet epc10bus,`
`uint16_t Command)`

Send command to the EPC10

Parameters

<i>epc10bus</i>	The EPC10 bus
<i>Command</i>	The command

10.37.3.13 SetADCOffset() void SetADCOffset (
 LIH30_ADC_Channel_EnumNet *AdcChannel*,
 uint32_t *Offset*)

Sets the ADC offset of the DACQ for a single channel

Parameters

<i>AdcChannel</i>	The ADC channel
<i>Offset</i>	The offset for the given channel number

10.37.3.14 SetDacIdleValue() void SetDacIdleValue (
 uint32_t *DacChannel*,
 int32_t *IdleValue*)

Sets the DAC Idle value

Parameters

<i>DacChannel</i>	The DAC channel
<i>IdleValue</i>	The idle value

10.37.3.15 SetDacUseIdleValue() void SetDacUseIdleValue (
 uint32_t *DacChannel*,
 bool *UseIdle*)

Sets if the DAC Idle value is used after stimulation

Parameters

<i>DacChannel</i>	The DAC channel
<i>UseIdle</i>	Use idle value

10.37.3.16 SetDigOutState() void SetDigOutState (
 uint16_t *DigOutState*)

Writes to the LIH30 digital output

Parameters

<i>DigOutState</i>	The bit mask defining the digital output state
--------------------	--

10.37.3.17 SetEepromPage() void SetEepromPage (
 uint32_t *EepromStartAddress*,
 array< int8_t >^ *EepromData*,
 LIH30_EPC10_Bus_EnumNet *epc10bus*)

Writes into EEprom on the EPC10 EEPROM

Parameters

<i>EepromStartAddress</i>	start address of memory area to write to
<i>EepromData</i>	pointer to internal memory for the supported amount of data
<i>epc10bus</i>	The EPC10 bus

10.37.3.18 SetSampleInterval() void SetSampleInterval (
 uint32_t *SampleInterval*)

Sets the Sample Interval for the DACQ and Stimulation

Parameters


<i>SampleInterval</i>	between the samples, Sample interval is available from 1 to 4194303
-----------------------	---

10.37.4 Property Documentation

10.37.4.1 StimulusFunction CStimulusFunctionNet^ StimulusFunction [get]

10.38 CMcsBus_AxisParametersNet Class Reference

Inheritance diagram for CMcsBus_AxisParametersNet:

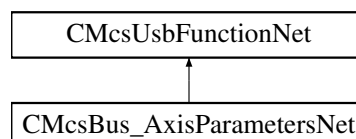

Public Member Functions

- [CMcsBus_AxisParametersNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_AxisParametersNet](#) (void)
- void [SetAxisParametersEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, unsigned int parameter)
- void [SetAxisParametersEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, int parameter)
- unsigned int [GetAxisParametersUnsignedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)
- int [GetAxisParametersSignedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)

Additional Inherited Members

10.38.1 Constructor & Destructor Documentation

10.38.1.1 [CMcsBus_AxisParametersNet\(\)](#) [CMcsBus_AxisParametersNet](#) (
[CMcsUsbNet](#)[^] device)

10.38.1.2 [~CMcsBus_AxisParametersNet\(\)](#) [~CMcsBus_AxisParametersNet](#) (
void)

10.38.2 Member Function Documentation

10.38.2.1 [GetAxisParametersSignedEeprom\(\)](#) int [GetAxisParametersSignedEeprom](#) (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned char *axis*,
unsigned short *index*)


10.38.2.2 [GetAxisParametersUnsignedEeprom\(\)](#) unsigned int [GetAxisParametersUnsignedEeprom](#) (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned char *axis*,
unsigned short *index*)

10.38.2.3 SetAxisParametersEeprom() [1/2] void SetAxisParametersEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *index*,
 int *parameter*)

10.38.2.4 SetAxisParametersEeprom() [2/2] void SetAxisParametersEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *index*,
 unsigned int *parameter*)

10.39 CMcsBus_ExtensionNet Class Reference

Inheritance diagram for CMcsBus_ExtensionNet:

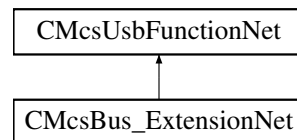

Public Member Functions

- [CMcsBus_ExtensionNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_ExtensionNet](#) (void)
- void [SetLEDSwitch](#) (unsigned char busnumber, unsigned char busaddress, unsigned short LEDSwitch)
- unsigned short [GetLEDSwitch](#) (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

10.39.1 Constructor & Destructor Documentation

10.39.1.1 CMcsBus_ExtensionNet() [CMcsBus_ExtensionNet](#) (
 [CMcsUsbNet](#)[^] device)

10.39.1.2 ~CMcsBus_ExtensionNet() [~CMcsBus_ExtensionNet](#) (
 void)


10.39.2 Member Function Documentation

10.39.2.1 GetLEDSwitch() unsigned short GetLEDSwitch (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.39.2.2 SetLEDSwitch() void SetLEDSwitch (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *LEDSwitch*)

10.40 CMcsBus_FYIExtensionNet Class Reference

Inheritance diagram for CMcsBus_FYIExtensionNet:

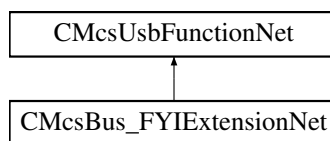

Public Member Functions

- [CMcsBus_FYIExtensionNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_FYIExtensionNet](#) (void)
- void [SetValves](#) (unsigned char busnumber, unsigned char busaddress, unsigned int states)
- unsigned int [GetValves](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetDIO](#) (unsigned char busnumber, unsigned char busaddress, unsigned short io)
- unsigned short [GetDIO](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetSingleHeater](#) (unsigned char busnumber, unsigned char busaddress, short index, unsigned short power)
- unsigned short [GetSingleHeater](#) (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members

10.40.1 Constructor & Destructor Documentation

10.40.1.1 CMcsBus_FYIExtensionNet() [CMcsBus_FYIExtensionNet](#) (
 [CMcsUsbNet](#)[^] *device*)

10.40.1.2 `~CMcsBus_FYIExtensionNet()` `~CMcsBus_FYIExtensionNet` (
 `void`)

10.40.2 Member Function Documentation

10.40.2.1 `GetDIO()` `unsigned short GetDIO` (
 `unsigned char busnumber`,
 `unsigned char busaddress`)

10.40.2.2 `GetSingleHeater()` `unsigned short GetSingleHeater` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `short index`)

10.40.2.3 `GetValves()` `unsigned int GetValves` (
 `unsigned char busnumber`,
 `unsigned char busaddress`)


10.40.2.4 `SetDIO()` `void SetDIO` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `unsigned short io`)

10.40.2.5 `SetSingleHeater()` `void SetSingleHeater` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `short index`,
 `unsigned short power`)

10.40.2.6 `SetValves()` `void SetValves` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `unsigned int states`)

10.41 CMcsBus_MotorControlNet Class Reference

Inheritance diagram for CMcsBus_MotorControlNet:

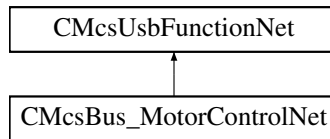

Public Member Functions

- [CMcsBus_MotorControlNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_MotorControlNet](#) (void)
- void [SetMCScalingFactorEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int [GetMCScalingFactorEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCScalingFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int [GetMCScalingFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCMaxSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCMaxSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravelEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravelEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCMaxCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCMaxCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRegulatorGainEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short [GetMCRegulatorGainEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRegulatorGain](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short [GetMCRegulatorGain](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCMaxAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void [SetMCMaxAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCMaxAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short [GetMCStandbyCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short [GetMCStandbyCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyTimeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short [GetMCStandbyTimeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyTime](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short [GetMCStandbyTime](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCBreakCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCBreakCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCBreakCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCBreakCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCConfigEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short [GetMCConfigEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCConfig](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short [GetMCConfig](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short [GetMCSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReferenceCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCReferenceCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReferenceCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCReferenceCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet [GetMCCurrentModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void [SetMCCurrentMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet [GetMCCurrentMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAxisRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short revision)
- unsigned short [GetMCAxisRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedUnitEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int32_t speedunit)
- int32_t [GetMCSpeedUnitEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, bool OnOff_status)
- bool [GetMCOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short [GetMCSpeedShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAccelerationShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAccelerationShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrentShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravelShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravelShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int [GetMCCurrentPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCNewPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int [GetMCNewPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- short [GetMCCurrentSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [StartMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRotation](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char onoff)
- unsigned short [GetMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReference](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char switch_enable, unsigned char switch_polarity)
- unsigned char [GetMCReference](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned char% switch_port)
- void [StopMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentModeShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet [GetMCCurrentModeShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short [GetMCPhase](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short [GetMCPhaseOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetSubChannel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short subchannel)
- unsigned short [GetSubChannel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

Additional Inherited Members

10.41.1 Constructor & Destructor Documentation

10.41.1.1 CMcsBus_MotorControlNet() `CMcsBus_MotorControlNet (CMcsUsbNet^ device)`

10.41.1.2 ~CMcsBus_MotorControlNet() `~CMcsBus_MotorControlNet (void)`

10.41.2 Member Function Documentation

10.41.2.1 GetMCAcceleration() `unsigned short GetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

10.41.2.2 GetMCAccelerationEeprom() `unsigned short GetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

10.41.2.3 GetMCAccelerationShortCommand() `unsigned short GetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

10.41.2.4 GetMCAxisRevisionEeprom() `unsigned short GetMCAxisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

10.41.2.5 GetMCBreakCurrent() short GetMCBreakCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.6 GetMCBreakCurrentEeprom() short GetMCBreakCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.7 GetMCConfig() unsigned short GetMCConfig (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.8 GetMCConfigEeprom() unsigned short GetMCConfigEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.9 GetMCCurrent() short GetMCCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.10 GetMCCurrentEeprom() short GetMCCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.11 GetMCCurrentMode() RoboCurrentModeEnumNet GetMCCurrentMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

- 10.41.2.12 GetMCCurrentModeEeprom()** RoboCurrentModeEnumNet GetMCCurrentModeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.13 GetMCCurrentModeShortCommand()** RoboCurrentModeEnumNet GetMCCurrentModeShort←
Command (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.14 GetMCCurrentPosition()** int GetMCCurrentPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.15 GetMCCurrentShortCommand()** short GetMCCurrentShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.16 GetMCCurrentSpeed()** short GetMCCurrentSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.17 GetMCMaxAcceleration()** unsigned short GetMCMaxAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.18 GetMCMaxAccelerationEeprom()** unsigned short GetMCMaxAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.19 GetMCMaxCurrent() short GetMCMaxCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.20 GetMCMaxCurrentEeprom() short GetMCMaxCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.21 GetMCMaxSpeed() unsigned short GetMCMaxSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.22 GetMCMaxSpeedEeprom() unsigned short GetMCMaxSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.23 GetMCMaxTravel() int GetMCMaxTravel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.24 GetMCMaxTravelEeprom() int GetMCMaxTravelEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.25 GetMCMaxTravelShortCommand() int GetMCMaxTravelShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

- 10.41.2.26 GetMCMovement()** unsigned short GetMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.27 GetMCNewPosition()** int GetMCNewPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.28 GetMCOutputOnOff()** bool GetMCOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.29 GetMCPhase()** unsigned short GetMCPhase (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.30 GetMCPhaseOffset()** unsigned short GetMCPhaseOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 10.41.2.31 GetMCReference()** unsigned char GetMCReference (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned char% *switch_port*)
- 10.41.2.32 GetMCReferenceCurrent()** short GetMCReferenceCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.33 GetMCReferenceCurrentEeprom() short GetMCReferenceCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.34 GetMCRegulatorGain() short GetMCRegulatorGain (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.35 GetMCRegulatorGainEeprom() short GetMCRegulatorGainEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.36 GetMCScalingFactor() int GetMCScalingFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.37 GetMCScalingFactorEeprom() int GetMCScalingFactorEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.38 GetMCSpeed() short GetMCSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.39 GetMCSpeedEeprom() unsigned short GetMCSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.40 GetMCSpeedShortCommand() short GetMCSpeedShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.41 GetMCSpeedUnitEeprom() int32_t GetMCSpeedUnitEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.42 GetMCStandbyCurrent() short GetMCStandbyCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.43 GetMCStandbyCurrentEeprom() short GetMCStandbyCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.44 GetMCStandbyTime() short GetMCStandbyTime (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.45 GetMCStandbyTimeEeprom() short GetMCStandbyTimeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.46 GetSubChannel() unsigned short GetSubChannel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.47 SetMCAcceleration() void SetMCAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

10.41.2.48 SetMCAccelerationEeprom() void SetMCAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

10.41.2.49 SetMCAccelerationShortCommand() void SetMCAccelerationShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

10.41.2.50 SetMCAxisRevisionEeprom() void SetMCAxisRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *revision*)

10.41.2.51 SetMCBreakCurrent() void SetMCBreakCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.52 SetMCBreakCurrentEeprom() void SetMCBreakCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.53 SetMCConfig() void SetMCConfig (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *config*)

10.41.2.54 SetMCConfigEeprom() void SetMCConfigEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *config*)

10.41.2.55 SetMCCurrent() void SetMCCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.56 SetMCCurrentEeprom() void SetMCCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.57 SetMCCurrentMode() void SetMCCurrentMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

10.41.2.58 SetMCCurrentModeEeprom() void SetMCCurrentModeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

10.41.2.59 SetMCCurrentModeShortCommand() void SetMCCurrentModeShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

10.41.2.60 SetMCCurrentPosition() void SetMCCurrentPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *position*)

10.41.2.61 SetMCCurrentShortCommand() void SetMCCurrentShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.62 SetMCMaxAcceleration() void SetMCMaxAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

10.41.2.63 SetMCMaxAccelerationEeprom() void SetMCMaxAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

10.41.2.64 SetMCMaxCurrent() void SetMCMaxCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.65 SetMCMaxCurrentEeprom() void SetMCMaxCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.66 SetMCMaxSpeed() void SetMCMaxSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

10.41.2.67 SetMCMaxSpeedEeprom() void SetMCMaxSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

10.41.2.68 SetMCMaxTravel() void SetMCMaxTravel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

10.41.2.69 SetMCMaxTravelEeprom() void SetMCMaxTravelEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

10.41.2.70 SetMCMaxTravelShortCommand() void SetMCMaxTravelShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

10.41.2.71 SetMCNewPosition() void SetMCNewPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *position*)

10.41.2.72 SetMCOutputOnOff() void SetMCOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 bool *OnOff_status*)

10.41.2.73 SetMCReference() void SetMCReference (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned char *switch_enable*,
 unsigned char *switch_polarity*)

10.41.2.74 SetMCReferenceCurrent() void SetMCReferenceCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.75 SetMCReferenceCurrentEeprom() void SetMCReferenceCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

10.41.2.76 SetMCRegulatorGain() void SetMCRegulatorGain (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *gain*)

10.41.2.77 SetMCRegulatorGainEeprom() void SetMCRegulatorGainEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *gain*)

10.41.2.78 SetMCRotation() void SetMCRotation (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned char *onoff*)

10.41.2.79 SetMCScalingFactor() void SetMCScalingFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *factor*)

10.41.2.80 SetMCScalingFactorEeprom() void SetMCScalingFactorEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *factor*)

10.41.2.81 SetMCSpeed() void SetMCSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *speed*)

10.41.2.82 SetMCSpeedEeprom() void SetMCSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

10.41.2.83 SetMCSpeedShortCommand() void SetMCSpeedShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *speed*)

10.41.2.84 SetMCSpeedUnitEeprom() void SetMCSpeedUnitEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int32_t *speedunit*)

10.41.2.85 SetMCStandbyCurrent() void SetMCStandbyCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *percent*)

10.41.2.86 SetMCStandbyCurrentEeprom() void SetMCStandbyCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *percent*)

10.41.2.87 SetMCStandbyTime() void SetMCStandbyTime (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *t*)

10.41.2.88 SetMCStandbyTimeEeprom() void SetMCStandbyTimeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *t*)


10.41.2.89 SetSubChannel() void SetSubChannel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *subchannel*)

10.41.2.90 StartMCMovement() void StartMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.41.2.91 StopMCMovement() void StopMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

10.42 CMcsBus_SensorNet Class Reference

Inheritance diagram for CMcsBus_SensorNet:

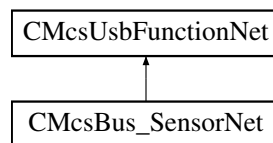

Public Member Functions

- [CMcsBus_SensorNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_SensorNet](#) (void)
- void [SetMinimalThreshold](#) (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short [GetMinimalThreshold](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetDetectionThreshold](#) (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short [GetDetectionThreshold](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetLatency](#) (unsigned char busnumber, unsigned char busaddress, unsigned short latency)
- unsigned short [GetLatency](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetBubbleStatus](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetLatencyCounter](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetDetectorValue](#) (unsigned char busnumber, unsigned char busaddress)
- array< int > [^] [GetPressure](#) (unsigned char busnumber, unsigned char busaddress, int n)
- int [GetPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRegulatorOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned char onoff)
- unsigned char [GetRegulatorOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)

- void [SetSollPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, int pressure)
- int [GetSollPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRegulatorFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, int factor)
- int [GetRegulatorFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- array< unsigned short > ^ [GetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress)
- int [GetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- unsigned int [GetRegulatorStatus](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRotatePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, short speed)
- short [GetRotatePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetMovePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned short speed, int position)
- void [GetMovePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% speed, [System::Runtime::InteropServices::Out]int% position)
- void [SetRegulationTimeouts](#) (unsigned char busnumber, unsigned char busaddress, unsigned short MaxSpeedWait, unsigned short MaxSignChange)
- void [GetRegulationTimeouts](#) (unsigned char busnumber, unsigned char busaddress, [System::Runtime::InteropServices::Out]unsigned short% MaxSpeedWait, [System::Runtime::InteropServices::Out]unsigned short% MaxSignChange)
- array< int > ^ [Get4ADC](#) (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ [Get4ADCAverage](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4DAC](#) (unsigned char busnumber, unsigned char busaddress, array< unsigned short > ^ dac)
- array< unsigned short > ^ [Get4DAC](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4ADCMode](#) (unsigned char busnumber, unsigned char busaddress, PatchServAdcModeEnumNet mode)
- PatchServAdcModeEnumNet [Get4ADCMode](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4ADCCatchampAverageShift](#) (unsigned char busnumber, unsigned char busaddress, unsigned int shift)
- unsigned int [Get4ADCCatchampAverageShift](#) (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ [Get2AnalogInput](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [Get2DigitalInput](#) (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ [GetADCs](#) (unsigned char busnumber, unsigned char busaddress, int n)
- array< unsigned short > ^ [GetADCsLoop](#) (unsigned char busnumber, unsigned char busaddress, int n)
- void [SetPiezoState](#) (unsigned char busnumber, unsigned char busaddress, int state)
- void [GetPiezoState](#) (unsigned char busnumber, unsigned char busaddress, [System::Runtime::InteropServices::Out]int% state, [System::Runtime::InteropServices::Out]int% reason)
- void [SetDACs](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, array< unsigned short > ^ dac_times_voltages)
- array< unsigned short > ^ [GetDACs](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetSamplePeriode](#) (unsigned char busnumber, unsigned char busaddress, unsigned short periode)
- unsigned short [GetSamplePeriode](#) (unsigned char busnumber, unsigned char busaddress)
- void [StartSync](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetSyncState](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacAmplitude](#) (unsigned char busnumber, unsigned char busaddress, unsigned short dacAmplitude)
- unsigned short [CatchAmpGetDacAmplitude](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacOffset](#) (unsigned char busnumber, unsigned char busaddress, short dacOffset)
- short [CatchAmpGetDacOffset](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcMean](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcValue](#) (unsigned char busnumber, unsigned char busaddress)

- int [CatchAmpGetAdcValueH](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcValueL](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetPwmEnable](#) (unsigned char busnumber, unsigned char busaddress, bool pwmEnable)
- bool [CatchAmpGetPwmEnable](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacEnable](#) (unsigned char busnumber, unsigned char busaddress, bool dacEnable)
- bool [CatchAmpGetDacEnable](#) (unsigned char busnumber, unsigned char busaddress)
- int [TactSwitchGetState](#) (unsigned char busnumber, unsigned char busaddress)
- void [TactSwitchSetDisplay](#) (unsigned char busnumber, unsigned char busaddress, int Melody)

Additional Inherited Members

10.42.1 Constructor & Destructor Documentation

10.42.1.1 CMcsBus_SensorNet() `CMcsBus_SensorNet (CMcsUsbNet^ device)`

10.42.1.2 ~CMcsBus_SensorNet() `~CMcsBus_SensorNet (void)`

10.42.2 Member Function Documentation

10.42.2.1 CatchAmpGetAdcMean() `int CatchAmpGetAdcMean (unsigned char busnumber, unsigned char busaddress)`

10.42.2.2 CatchAmpGetAdcValue() `int CatchAmpGetAdcValue (unsigned char busnumber, unsigned char busaddress)`

10.42.2.3 CatchAmpGetAdcValueH() `int CatchAmpGetAdcValueH (unsigned char busnumber, unsigned char busaddress)`

10.42.2.4 CatchAmpGetAdcValueL() int CatchAmpGetAdcValueL (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.5 CatchAmpGetDacAmplitude() unsigned short CatchAmpGetDacAmplitude (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.6 CatchAmpGetDacEnable() bool CatchAmpGetDacEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.7 CatchAmpGetDacOffset() short CatchAmpGetDacOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.8 CatchAmpGetPwmEnable() bool CatchAmpGetPwmEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.9 CatchAmpSetDacAmplitude() void CatchAmpSetDacAmplitude (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *dacAmplitude*)

10.42.2.10 CatchAmpSetDacEnable() void CatchAmpSetDacEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 bool *dacEnable*)

10.42.2.11 CatchAmpSetDacOffset() void CatchAmpSetDacOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *dacOffset*)

10.42.2.12 CatchAmpSetPwmEnable() void CatchAmpSetPwmEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 bool *pwmEnable*)

10.42.2.13 Get2AnalogInput() array<unsigned short> ^ Get2AnalogInput (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.14 Get2DigitalInput() unsigned short Get2DigitalInput (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.15 Get4ADC() array<int> ^ Get4ADC (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.16 Get4ADCAverage() array<int> ^ Get4ADCAverage (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.17 Get4ADCCatchampAverageShift() unsigned int Get4ADCCatchampAverageShift (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.18 Get4ADCMode() PatchServAdcModeEnumNet Get4ADCMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.19 Get4DAC() array<unsigned short> ^ Get4DAC (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.20 GetADCs() array<unsigned short> ^ GetADCs (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *n*)

10.42.2.21 GetADCsLoop() array<unsigned short> ^ GetADCsLoop (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *n*)

10.42.2.22 GetBubbleStatus() unsigned short GetBubbleStatus (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.23 GetDACs() array<unsigned short> ^ GetDACs (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.24 GetDetectionThreshold() unsigned short GetDetectionThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.25 GetDetectorValue() unsigned short GetDetectorValue (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.26 GetLatency() unsigned short GetLatency (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.27 GetLatencyCounter() unsigned short GetLatencyCounter (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.28 GetMinimalThreshold() unsigned short GetMinimalThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.29 GetMovePump() void GetMovePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned short% *speed*,
 [System::Runtime::InteropServices::Out] int% *position*)

10.42.2.30 GetPiezoState() void GetPiezoState (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 [System::Runtime::InteropServices::Out] int% *state*,
 [System::Runtime::InteropServices::Out] int% *reason*)

10.42.2.31 GetPressure() [1/2] array<int> ^ GetPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *n*)

10.42.2.32 GetPressure() [2/2] int GetPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.33 GetPressureOffset() [1/2] array<unsigned short> ^ GetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.34 GetPressureOffset() [2/2] int GetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.35 GetRegulationTimeouts() void GetRegulationTimeouts (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 [System::Runtime::InteropServices::Out] unsigned short% *MaxSpeedWait*,
 [System::Runtime::InteropServices::Out] unsigned short% *MaxSignChange*)

10.42.2.36 GetRegulatorFactor() int GetRegulatorFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.37 GetRegulatorOnOff() unsigned char GetRegulatorOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.38 GetRegulatorStatus() unsigned int GetRegulatorStatus (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.39 GetRotatePump() short GetRotatePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.40 GetSamplePeriode() unsigned short GetSamplePeriode (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.41 GetSollPressure() int GetSollPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.42 GetSyncState() unsigned short GetSyncState (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.43 Set4ADCCatchampAverageShift() void Set4ADCCatchampAverageShift (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned int *shift*)

10.42.2.44 Set4ADCMode() void Set4ADCMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 PatchServAdcModeEnumNet *mode*)

10.42.2.45 Set4DAC() void Set4DAC (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 array< unsigned short >^ *dac*)

10.42.2.46 SetDACs() void SetDACs (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 array< unsigned short >^ *dac_times_voltages*)

10.42.2.47 SetDetectionThreshold() void SetDetectionThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *threshold*)

10.42.2.48 SetLatency() void SetLatency (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *latency*)

10.42.2.49 SetMinimalThreshold() void SetMinimalThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *threshold*)

10.42.2.50 SetMovePump() void SetMovePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 unsigned short *speed*,
 int *position*)

10.42.2.51 SetPiezoState() void SetPiezoState (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *state*)

10.42.2.52 SetPressureOffset() void SetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

10.42.2.53 SetRegulationTimeouts() void SetRegulationTimeouts (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *MaxSpeedWait*,
 unsigned short *MaxSignChange*)

10.42.2.54 SetRegulatorFactor() void SetRegulatorFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 int *factor*)

10.42.2.55 SetRegulatorOnOff() void SetRegulatorOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 unsigned char *onoff*)

10.42.2.56 SetRotatePump() void SetRotatePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 short *speed*)

10.42.2.57 SetSamplePeriode() void SetSamplePeriode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *periode*)

10.42.2.58 SetSollPressure() void SetSollPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 int *pressure*)


10.42.2.59 StartSync() void StartSync (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.60 TactSwitchGetState() int TactSwitchGetState (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.42.2.61 TactSwitchSetDisplay() void TactSwitchSetDisplay (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *Melody*)

10.43 CMcsBus_TempSensorNet Class Reference

Inheritance diagram for CMcsBus_TempSensorNet:

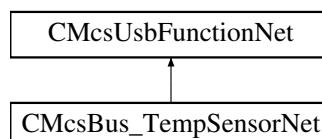

Public Member Functions

- [CMcsBus_TempSensorNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_TempSensorNet](#) (void)
- short [GetTemperatur](#) (unsigned char busnumber, unsigned char busaddress)
- short [GetTemperatur](#) (unsigned char busnumber, unsigned char busaddress, short index)
- void [SetNanoVoltsPerKelvin](#) (unsigned char busnumber, unsigned char busaddress, int nanovoltsperkelvin)
- int [GetNanoVoltsPerKelvin](#) (unsigned char busnumber, unsigned char busaddress)
- short [GetThermoVoltage](#) (unsigned char busnumber, unsigned char busaddress, short index)
- short [GetThermoTemp](#) (unsigned char busnumber, unsigned char busaddress, short index)
- void [SetThermoOffset](#) (unsigned char busnumber, unsigned char busaddress, short index, short offset)
- short [GetThermoOffset](#) (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members**10.43.1 Constructor & Destructor Documentation**

10.43.1.1 [CMcsBus_TempSensorNet\(\)](#) [CMcsBus_TempSensorNet](#) (
[CMcsUsbNet](#)[^] device)

10.43.1.2 [~CMcsBus_TempSensorNet\(\)](#) [~CMcsBus_TempSensorNet](#) (
void)

10.43.2 Member Function Documentation

10.43.2.1 [GetNanoVoltsPerKelvin\(\)](#) int [GetNanoVoltsPerKelvin](#) (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.43.2.2 [GetTemperatur\(\)](#) [1/2] short [GetTemperatur](#) (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.43.2.3 [GetTemperatur\(\)](#) [2/2] short [GetTemperatur](#) (
unsigned char *busnumber*,
unsigned char *busaddress*,
short *index*)

10.43.2.4 GetThermoOffset() `short GetThermoOffset (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `short index)`

10.43.2.5 GetThermoTemp() `short GetThermoTemp (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `short index)`


10.43.2.6 GetThermoVoltage() `short GetThermoVoltage (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `short index)`

10.43.2.7 SetNanoVoltsPerKelvin() `void SetNanoVoltsPerKelvin (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `int nanovoltsperkelvin)`

10.43.2.8 SetThermoOffset() `void SetThermoOffset (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `short index,`
 `short offset)`

10.44 CMcsBus_VoltageModeNet Class Reference

Inheritance diagram for CMcsBus_VoltageModeNet:

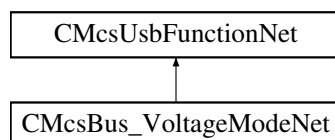

Public Member Functions

- [CMcsBus_VoltageModeNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_VoltageModeNet](#) (void)
- void [SetVMMMaxPositiveCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxPositiveCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxPositiveCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxNegativeCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxNegativeCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxPositiveVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxPositiveVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxNegativeVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxNegativeVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMOOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, unsigned short status)
- unsigned short [GetVMOOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)

Additional Inherited Members**10.44.1 Constructor & Destructor Documentation**
10.44.1.1 CMcsBus_VoltageModeNet() [CMcsBus_VoltageModeNet](#) ([CMcsUsbNet](#)[^] device)

10.44.1.2 `~CMcsBus_VoltageModeNet()` `~CMcsBus_VoltageModeNet` (
 void)

10.44.2 Member Function Documentation

10.44.2.1 `GetVMMaxNegativeCurrent()` short `GetVMMaxNegativeCurrent` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.2 `GetVMMaxNegativeCurrentEeprom()` short `GetVMMaxNegativeCurrentEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.3 `GetVMMaxNegativeVoltage()` short `GetVMMaxNegativeVoltage` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.4 `GetVMMaxNegativeVoltageEeprom()` short `GetVMMaxNegativeVoltageEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.5 `GetVMMaxPositiveCurrent()` short `GetVMMaxPositiveCurrent` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.6 `GetVMMaxPositiveCurrentEeprom()` short `GetVMMaxPositiveCurrentEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.7 GetVMMaxPositiveVoltage() short GetVMMaxPositiveVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.8 GetVMMaxPositiveVoltageEeprom() short GetVMMaxPositiveVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.9 GetVMOutputOnOff() unsigned short GetVMOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.10 GetVMVoltage() short GetVMVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

10.44.2.11 SetVMMaxNegativeCurrent() void SetVMMaxNegativeCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

10.44.2.12 SetVMMaxNegativeCurrentEeprom() void SetVMMaxNegativeCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

10.44.2.13 SetVMMaxNegativeVoltage() void SetVMMaxNegativeVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

10.44.2.14 SetVMMaxNegativeVoltageEeprom() void SetVMMaxNegativeVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

10.44.2.15 SetVMMaxPositiveCurrent() void SetVMMaxPositiveCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

10.44.2.16 SetVMMaxPositiveCurrentEeprom() void SetVMMaxPositiveCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

10.44.2.17 SetVMMaxPositiveVoltage() void SetVMMaxPositiveVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)


10.44.2.18 SetVMMaxPositiveVoltageEeprom() void SetVMMaxPositiveVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

10.44.2.19 SetVMOutputOnOff() void SetVMOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 unsigned short *status*)

10.44.2.20 SetVMVoltage() void SetVMVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

10.45 CMcsBusNet Class Reference

Inheritance diagram for CMcsBusNet:

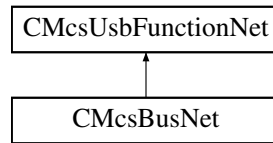

Public Member Functions

- [CMcsBusNet](#) ([CMcsUsbNet](#)[^] device)
- virtual [~CMcsBusNet](#) (void)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, short value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned int value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, int value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned short% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]short% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned int% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]int% value)
- void [SetBusAddressEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short [GetBusAddressEeprom](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetBusAddress](#) (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short [GetBusAddress](#) (unsigned char busnumber, unsigned char busaddress)
- void [CMcsBusNet::SetModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short [CMcsBusNet::GetModeEeprom](#) (unsigned char busnumber, unsigned char busaddress)
- void [CMcsBusNet::SetMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short [CMcsBusNet::GetMode](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetHWRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short revision)
- unsigned short [GetHWRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

10.45.1 Constructor & Destructor Documentation

10.45.1.1 CMcsBusNet() [CMcsBusNet](#) ([CMcsUsbNet](#)[^] device)

10.45.1.2 `~CMcsBusNet()` virtual `~CMcsBusNet` (
void) [virtual]

10.45.2 Member Function Documentation

10.45.2.1 `CMcsBusNet::GetMode()` unsigned short `CMcsBusNet::GetMode` (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.45.2.2 `CMcsBusNet::GetModeEeprom()` unsigned short `CMcsBusNet::GetModeEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.45.2.3 `CMcsBusNet::SetMode()` void `CMcsBusNet::SetMode` (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned short *mode*)

10.45.2.4 `CMcsBusNet::SetModeEeprom()` void `CMcsBusNet::SetModeEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned short *mode*)

10.45.2.5 `GetBusAddress()` unsigned short `GetBusAddress` (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.45.2.6 `GetBusAddressEeprom()` unsigned short `GetBusAddressEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*)

10.45.2.7 GetCommand() [1/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] int% *value*)

10.45.2.8 GetCommand() [2/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] short% *value*)

10.45.2.9 GetCommand() [3/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned int% *value*)

10.45.2.10 GetCommand() [4/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned short% *value*)

10.45.2.11 GetHWRevisionEeprom() unsigned short GetHWRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

10.45.2.12 SetBusAddress() void SetBusAddress (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *newaddress*)

10.45.2.13 SetBusAddressEeprom() void SetBusAddressEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *newaddress*)

10.45.2.14 SetCommand() [1/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *value*)

10.45.2.15 SetCommand() [2/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *value*)

10.45.2.16 SetCommand() [3/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned int *value*)


10.45.2.17 SetCommand() [4/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *value*)

10.45.2.18 SetHWRevisionEeprom() void SetHWRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *revision*)

10.46 CMcsUsbDacqNet Class Reference

Base class for data acquisition devices.

Inheritance diagram for CMcsUsbDacqNet:

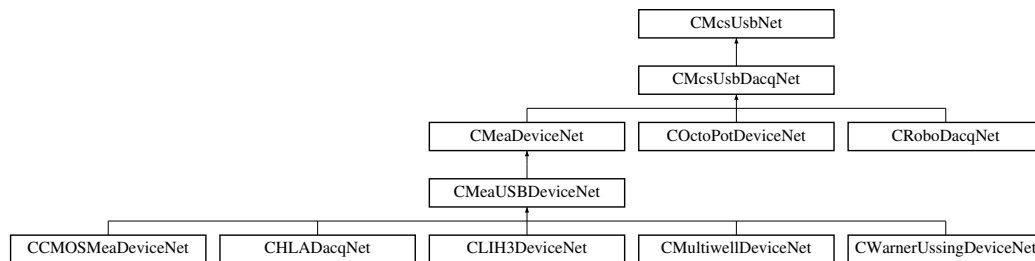

Classes

- class [CHWInfo](#)
Class to provide hardware information about the device.

Public Member Functions

- [CMcsUsbDacqNet](#) ()
- [~CMcsUsbDacqNet](#) ()
- virtual uint32_t [GetVoltageRangeIndex](#) (unsigned int virtualDevice)
- virtual void [SetVoltageRangeByIndex](#) (int32_t voltageRangeIndex, unsigned int virtualDevice)
Sets the voltage range on devices which support multiple voltage ranges.
- virtual void [SetVoltageRangeInMicroVolt](#) (int32_t voltageRange, unsigned int virtualDevice)
Sets the voltage range on devices which support multiple voltage ranges.
- virtual int32_t [GetVoltageRangeInMicroVolt](#) (unsigned int virtualDevice)
Gets the currently selected voltage range on devices which support multiple voltage ranges.
- virtual int32_t [GetVoltageRangeInMilliVolt](#) ()
Gets the currently selected voltage range on devices which support multiple voltage ranges.
- virtual void [SetDataMode](#) (DataModeEnumNet dataMode, unsigned int virtualDevice)
Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.
- virtual DataModeEnumNet [GetDataMode](#) (unsigned int virtualDevice)
Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.
- void [SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- void [SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, W2100DigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- void [SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, SCUDigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- void [SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MEA2100_256DigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.

- template<typename digitalsourceenum >
void [SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [DigitalSource](#)< digitalsourceenum >^ source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- void [GetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
Gets the function/source of an digital output bit.
- void [GetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
Gets the function/source of an digital output bit.
- void [GetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]SCUDigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
Gets the function/source of an digital output bit.
- void [GetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]MEA2100_256DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
Gets the function/source of an digital output bit.
- template<typename digitalsourceenum >
void [GetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out][DigitalSource](#)< digitalsourceenum >^% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
Gets the function/source of an digital output bit.
- virtual AdapterTypeEnumNet [GetAdapterType](#) ()
Gets the adapter which is connected to the MEA2100 device.
- virtual MeaLayoutEnumNet [GetMeaLayout](#) ()
Gets the MEA layout which is connected to the MEA2100 device.
- virtual uint32_t [GetAdcDataFormat](#) (uint32_t virtualDevice)
Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.
- virtual uint32_t [GetAnalogValueUnit](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] AnalogUnitEnumNet% unit)
- virtual uint32_t [GetResolutionPerDigit](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% res, [System::Runtime::InteropServices::Out] int% resUnit)
- virtual uint32_t [GetAdcZero](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% adcz)
- virtual uint32_t [GetHardwareMinRange](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t [GetHardwareMaxRange](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t [GetDataFormat](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)
- virtual uint32_t [GetNumberOfDataBits](#) (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)
Get the real number of data bits.
- virtual void [SetSamplerate](#) (int32_t rate, unsigned int oversample, unsigned int virtualDevice)
Sets the sampling frequency of the device.
- virtual int32_t [GetSamplerate](#) (unsigned int virtualDevice)
Gets the sampling frequency of the device.
- virtual uint32_t [GetMaxSamplingFrequency](#) (int virtualDevice)
Gets the maximal sampling frequency of the device.
- virtual uint32_t [GetMinSamplingFrequencyStepsize](#) ()
Gets the minimal sampling frequency step size increment value of the device.

- virtual int32_t [GetChannelsInBlock](#) (unsigned int virtualDevice)
Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.
- virtual uint32_t [GetChannelLayout](#) ([System::Runtime::InteropServices::Out]int% AnalogChannels, [System::Runtime::InteropServices::Out]int% DigitalChannels, [System::Runtime::InteropServices::Out]int% ChecksumChannels, [System::Runtime::InteropServices::Out]int% TimestampChannels, [System::Runtime::InteropServices::Out]int% ChannelsInBlock, unsigned int virtualDevice)
- virtual void [SendStartDacq](#) ()
Start sampling.
- virtual void [SendStartDacq](#) (int VirtualDacqMap)
Start sampling.
- virtual void [SendStartStgAndDacq](#) (uint32_t trigger_map, int VirtualDacqMap)
Start sampling together with the STG.
- virtual void [SendStopDacq](#) ()
Stop sampling.
- virtual void [SendStopDacq](#) (int VirtualDacqMap)
Stop sampling.

Parameters

VirtualDacqMap	
----------------	--

- virtual void [SendStopStgAndDacq](#) (uint32_t trigger_map, int VirtualDacqMap)
Stop sampling together with the STG.
- virtual void [SendStopStgAndDacqWithOptions](#) (uint32_t trigger_map, int VirtualDacqMap, int options)
Stop sampling together with the STG and options.
- virtual void [StartLoop](#) ()
Start the data acquisition thread.
- virtual void [StartLoop](#) (int32_t timeout)
Start the data acquisition thread.
- virtual void [StartLoop](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb)
Start the data acquisition thread.
- virtual void [StartLoop](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)
Start the data acquisition thread.
- virtual void [StopLoop](#) ()
- virtual void [ClearBuffers](#) ()
- virtual void [StartDacq](#) ()
Start the data acquisition thread and sampling.
- virtual void [StartDacq](#) (int32_t timeout)
Start the data acquisition thread and sampling.
- virtual void [StartDacq](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb)
Start the data acquisition thread and sampling.
- virtual void [StartDacq](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)
Start the data acquisition thread and sampling.
- virtual void [StopDacq](#) ()
Stop the data acquisition thread and sampling.
- virtual void [StopDacq](#) (uint32_t virtualDevice)
Stop the data acquisition thread and sampling.

- virtual uint32_t [SetPoti](#) (uint32_t channel, uint32_t value, bool write_nvram)
- virtual uint32_t [GetPoti](#) (uint32_t channel, [System::Runtime::InteropServices::Out]uint32_t% value)
- virtual [CFilterPropertyNet](#) ^ [GetFilterProperty](#) (DacqGroupChannelEnumNet GroupID, unsigned int index)
- virtual array< [CFilterPropertyNet](#) ^> ^ [CMcsUsbDacqNet::GetFilterProperties](#) (DacqGroupChannelEnumNet GroupID)
- int [GetChannelDataFillSize](#) ()
- virtual void [SetSelectedChannels](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void [SetSelectedChannels](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedChannels](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void [SetSelectedChannels](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedData](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void [SetSelectedData](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedData](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void [SetSelectedData](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize)
- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize)
- virtual void [SetSelectedChannelsQueue](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual void [SetSelectedChannelsQueue](#) (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedChannelsQueue](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual void [SetSelectedChannelsQueue](#) (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual uint32_t [ChannelBlock_AvailFrames](#) (int handle)

Get the number of sample frames already available in the FIFO.

- virtual uint32_t [ChannelBlock_AvailFrames](#) (int handle, int queue)
- virtual array< uint16_t > ^ [ChannelBlock_ReadFramesUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format

- virtual uint32_t [ChannelBlock_ReadFramesUI16](#) (int handle, array< uint16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format

- virtual array< int16_t > ^ [ChannelBlock_ReadFramesI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format

- virtual uint32_t [ChannelBlock_ReadFramesI16](#) (int handle, array< int16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format

- virtual array< uint32_t > ^ [ChannelBlock_ReadFramesUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

- virtual uint32_t [ChannelBlock_ReadFramesUI32](#) (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

- virtual array< int32_t > ^ [ChannelBlock_ReadFramesI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

- virtual uint32_t [ChannelBlock_ReadFramesI32](#) (int handle, array< int32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

- virtual array< array< uint16_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< uint16_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI16](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< int16_t >^> ^ [ChannelBlock_ReadAsFrameArrayI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< int16_t >^> ^ [ChannelBlock_ReadAsFrameArrayI16](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< uint32_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< uint32_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI32](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< int32_t >^> ^ [ChannelBlock_ReadAsFrameArrayI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual array< array< int32_t >^> ^ [ChannelBlock_ReadAsFrameArrayI32](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

- virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ [ChannelBlock_ReadFramesDictUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

- virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ [ChannelBlock_ReadFramesDictI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

- virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ [ChannelBlock_ReadFramesDictUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ [ChannelBlock_ReadFramesDictI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ [GetGroupChannelDataUI16](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ [GetGroupChannelDataI16](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ [GetGroupChannelDataUI32](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ [GetGroupChannelDataI32](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- void [SetupGroupDacqQueue](#) (int queuesize, int threshold)
- [CHWInfo](#) ^ [HWInfo](#) ()

Static Public Attributes

- static const int [Error_Callback_Queue_Full](#) = 0x100
- static const int [Error_Callback_Aquisition_Stopped](#) = 0x200
- static const int [Error_Callback_Packet_Error](#) = 1
- static const int [Error_Callback_RingQueue_Full](#) = 3
- static const int [Error_Callback_Frames_Lost](#) = 4
- static const int [Error_Callback_Data_lost](#) = 5

Properties

- virtual int [Samplerate](#) [get, set]
The sampling frequency of the device in Hz.

Events

- [OnChannelData](#) ^ [ChannelDataEvent](#)
- [OnError](#) ^ [ErrorEvent](#)

Additional Inherited Members

10.46.1 Detailed Description

Base class for data acquisition devices.

10.46.2 Constructor & Destructor Documentation

10.46.2.1 CMcsUsbDacqNet() `CMcsUsbDacqNet ()`

10.46.2.2 ~CMcsUsbDacqNet() `~CMcsUsbDacqNet ()`

10.46.3 Member Function Documentation

10.46.3.1 AddSelectedChannelsQueue() [1/4] `virtual int AddSelectedChannelsQueue (`
`int nByteOffset,`
`int nChannelOffset,`
`array< bool >^ selectedChannels,`
`int queuesize,`
`int threshold,`
`SampleSizeNet samplesize) [virtual]`

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use `ChannelBlock_Read↔FramesDict...` with `handle = 0` to read the data.

When using 32 bit data format, `ChannelsInBlock` is still the number of 16 bit channels per frame, as obtained from `GetChannelsInBlock`, while `nChannels` is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format `nChannels = ChannelsInBlock/2`

Parameters

<code>nByteOffset</code>	Number of bytes to start with.
--------------------------	--------------------------------

Parameters

<code>nChannelOffset</code>	Number of channel to start with (counted in <code>samplesize</code> bytes).
-----------------------------	---

Parameters

<code>selectedChannels</code>	List of channels to be collected in the FIFO.
-------------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Returns

The handle to the Queue.

10.46.3.2 AddSelectedChannelsQueue() [2/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*) [virtual]

10.46.3.3 AddSelectedChannelsQueue() [3/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*) [virtual]

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read↔FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>nByteOffset</i>	Number of bytes to start with.
--------------------	--------------------------------

Parameters

<i>nChannelOffset</i>	Number of channel to start with (counted in samplesize bytes).
-----------------------	--

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Returns

The handle to the Queue.

10.46.3.4 AddSelectedChannelsQueue() [4/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*) [virtual]

10.46.3.5 ChannelBlock_AvailFrames() [1/2] virtual uint32_t ChannelBlock_AvailFrames (
 int *handle*) [virtual]

Get the number of sample frames already available in the FIFO.

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Returns

Number of sample frames available in the FIFO.

10.46.3.6 ChannelBlock_AvailFrames() [2/2] virtual uint32_t ChannelBlock_AvailFrames (
 int *handle*,
 int *queue*) [virtual]

10.46.3.7 ChannelBlock_ReadAsFrameArrayI16() [1/2] virtual array<array<int16_t>> ^ Channel↔
 Block_ReadAsFrameArrayI16 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int16_t frame arrays.

10.46.3.8 ChannelBlock_ReadAsFrameArrayI16() [2/2] virtual array<array<int16_t>^> ^ Channel↔

```
Block_ReadAsFrameArrayI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int16_t frame arrays.

10.46.3.9 ChannelBlock_ReadAsFrameArrayI32() [1/2] virtual array<array<int32_t>^> ^ Channel↔
 Block_ReadAsFrameArrayI32 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int32_t frame arrays.

10.46.3.10 ChannelBlock_ReadAsFrameArrayI32() [2/2] virtual array<array<int32_t>^> ^ Channel↔
 Block_ReadAsFrameArrayI32 (
 int *handle*,
 int *queue*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int32_t frame arrays.

10.46.3.11 ChannelBlock_ReadAsFrameArrayUI16() [1/2] virtual array<array<uint16_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint16_t frame arrays.

10.46.3.12 ChannelBlock_ReadAsFrameArrayUI16() [2/2] virtual array<array<uint16_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint16_t frame arrays.

10.46.3.13 ChannelBlock_ReadAsFrameArrayUI32() [1/2] virtual array<array<uint32_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint32_t frame arrays.

10.46.3.14 ChannelBlock_ReadAsFrameArrayUI32() [2/2] virtual array<array<uint32_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI32 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint32_t frame arrays.

10.46.3.15 ChannelBlock_ReadFramesDictI16() virtual System::Collections::Generic::Dictionary<int, array<int16_t>^> ^ ChannelBlock_ReadFramesDictI16 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int16_t arrays and hardware channel as key.

10.46.3.16 ChannelBlock_ReadFramesDictI32() virtual System::Collections::Generic::Dictionary<int, array<int32_t>^> ^ ChannelBlock_ReadFramesDictI32 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int32_t arrays and hardware channel as key.

10.46.3.17 ChannelBlock_ReadFramesDictUI16() `virtual System::Collections::Generic::Dictionary<int, array<uint16_t>> ^ ChannelBlock_ReadFramesDictUI16 (int handle, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint16_t arrays and hardware channel as key.

10.46.3.18 ChannelBlock_ReadFramesDictUI32() virtual System::Collections::Generic::Dictionary<int, array<uint32_t>^> ^ ChannelBlock_ReadFramesDictUI32 (int *handle*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint32_t arrays and hardware channel as key.

10.46.3.19 ChannelBlock_ReadFramesI16() [1/2] virtual uint32_t ChannelBlock_ReadFramesI16 (int *handle*, array< int16_t >^ *buffer*, int *frames_pos*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in int16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Error Status. 0 on success.

10.46.3.20 ChannelBlock_ReadFramesI16() [2/2] `virtual array<int16_t> ^ ChannelBlock_ReadFramesI16 (`
`int handle,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in int16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

10.46.3.21 ChannelBlock_ReadFramesI32() [1/2] `virtual uint32_t ChannelBlock_ReadFramesI32 (`
`int handle,`
`array< int32_t >^ buffer,`
`int frames_pos,`

```
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Error Status. 0 on success.

10.46.3.22 ChannelBlock_ReadFramesI32() [2/2] virtual array<uint32_t> ^ ChannelBlock_ReadFramesI32 (

```
int handle,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

10.46.3.23 ChannelBlock_ReadFramesUI16() [1/2] `virtual uint32_t ChannelBlock_ReadFramesUI16 (int handle, array< uint16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Error Status. 0 on success.

10.46.3.24 ChannelBlock_ReadFramesUI16() [2/2] `virtual array<uint16_t> ^ ChannelBlock_ReadFramesUI16 (int handle, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of data from the device.

10.46.3.25 ChannelBlock_ReadFramesUI32() [1/2] `virtual uint32_t ChannelBlock_ReadFramesUI32 (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Error Status. 0 on success.

10.46.3.26 ChannelBlock_ReadFramesUI32() [2/2] `virtual array<uint32_t> ^ ChannelBlock_ReadFramesUI32 (`
`int handle,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

10.46.3.27 ClearBuffers() `virtual void ClearBuffers () [virtual]`

10.46.3.28 CMcsUsbDacqNet::GetFilterProperties() virtual array<CFilterPropertyNet^> ^ CMcsUsbDacqNet::GetFilterProperties (DacqGroupChannelEnumNet GroupID) [virtual]

10.46.3.29 GetAdapterType() virtual AdapterTypeEnumNet GetAdapterType () [virtual]

Gets the adapter which is connected to the MEA2100 device.

Returns

AdapterTypeEnumNet which enumerates the possible adapters.

10.46.3.30 GetAdcDataFormat() virtual uint32_t GetAdcDataFormat (uint32_t virtualDevice) [virtual]

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

Returns

The data format in bits.

10.46.3.31 GetAdcZero() virtual uint32_t GetAdcZero (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% adcz) [virtual]

10.46.3.32 GetAnalogValueUnit() virtual uint32_t GetAnalogValueUnit (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] AnalogUnitEnumNet% unit) [virtual]

10.46.3.33 GetChannelDataFillSize() int GetChannelDataFillSize ()

10.46.3.34 GetChannelLayout() virtual uint32_t GetChannelLayout (
[System::Runtime::InteropServices::Out] int% *AnalogChannels*,
[System::Runtime::InteropServices::Out] int% *DigitalChannels*,
[System::Runtime::InteropServices::Out] int% *ChecksumChannels*,
[System::Runtime::InteropServices::Out] int% *TimestampChannels*,
[System::Runtime::InteropServices::Out] int% *ChannelsInBlock*,
unsigned int *virtualDevice*) [virtual]

10.46.3.35 GetChannelsInBlock() virtual int32_t GetChannelsInBlock (
unsigned int *virtualDevice*) [virtual]

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

Returns

Number of 16 bit datawords per sample frame.

10.46.3.36 GetDataFormat() virtual uint32_t GetDataFormat (
uint32_t *virtualDevice*,
DacqGroupChannelEnumNet *group*,
[System::Runtime::InteropServices::Out] int% *numberOfBits*) [virtual]

10.46.3.37 GetDataMode() virtual DataModeEnumNet GetDataMode (
unsigned int *virtualDevice*) [virtual]

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------

Returns

DataModeEnumNet which enumerates the possible data modes.

10.46.3.38 GetDigitalSource() [1/5] void GetDigitalSource (
DigitalTargetEnumNet *digitaltarget*,
int32_t *NrChannel*,
[System::Runtime::InteropServices::Out] [DigitalSource](#)< digitalsourceenum >^%
source,
[System::Runtime::InteropServices::Out] int% *bitnumber_offset*)

Gets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.39 GetDigitalSource() [2/5] `void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out] DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out] int% bitnumber_offset)`

Gets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.40 GetDigitalSource() [3/5] `void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out] MEA2100_256DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out] int% bitnumber_offset)`

Gets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.41 GetDigitalSource() [4/5] `void GetDigitalSource (DigitalTargetEnumNet digitaltarget,`

```
int32_t NrChannel,
[System::Runtime::InteropServices::Out] SCUDigitalSourceEnumNet% source,
[System::Runtime::InteropServices::Out] int% bitnumber_offset )
```

Gets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.42 GetDigitalSource() [5/5] void GetDigitalSource (DigitalTargetEnumNet *digitaltarget*, int32_t *NrChannel*, [System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% *source*, [System::Runtime::InteropServices::Out] int% *bitnumber_offset*)

Gets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.43 GetFilterProperty() virtual CFilterPropertyNet ^ GetFilterProperty (DacqGroupChannelEnumNet *GroupID*, unsigned int *index*) [virtual]

10.46.3.44 GetGroupChannelDataI16() virtual System::Collections::Generic::Dictionary<int, array<int16_t>^> ^ GetGroupChannelDataI16 (DacqGroupChannelEnumNet *group*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int16_t arrays and hardware channel as key.

10.46.3.45 GetGroupChannelDataI32() `virtual System::Collections::Generic::Dictionary<int, array<int32_t>^> ^ GetGroupChannelDataI32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int32_t arrays and hardware channel as key.

10.46.3.46 GetGroupChannelDataUI16() virtual System::Collections::Generic::Dictionary<int, array<uint16_t>^> ^ GetGroupChannelDataUI16 (DacqGroupChannelEnumNet *group*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint16_t arrays and hardware channel as key.

10.46.3.47 GetGroupChannelDataUI32() virtual System::Collections::Generic::Dictionary<int, array<uint32_t>^> ^ GetGroupChannelDataUI32 (DacqGroupChannelEnumNet *group*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint32_t arrays and hardware channel as key.

10.46.3.48 GetHardwareMaxRange() virtual uint32_t GetHardwareMaxRange (
uint32_t *virtualDevice*,
DacqGroupChannelEnumNet *group*,
[System::Runtime::InteropServices::Out] int% *r*,
[System::Runtime::InteropServices::Out] int% *rUnit*) [virtual]

10.46.3.49 GetHardwareMinRange() virtual uint32_t GetHardwareMinRange (
uint32_t *virtualDevice*,
DacqGroupChannelEnumNet *group*,
[System::Runtime::InteropServices::Out] int% *r*,
[System::Runtime::InteropServices::Out] int% *rUnit*) [virtual]

10.46.3.50 GetMaxSamplingFrequency() virtual uint32_t GetMaxSamplingFrequency (
int *virtualDevice*) [virtual]

Gets the maximal sampling frequency of the device.

Returns

Sampling frequency in Hz.

10.46.3.51 GetMeaLayout() virtual MeaLayoutEnumNet GetMeaLayout () [virtual]

Gets the MEA layout which is connected to the MEA2100 device.

Returns

MeaLayoutEnumNet which enumerates the MEA types.

10.46.3.52 GetMinSamplingFrequencyStepsize() virtual uint32_t GetMinSamplingFrequencyStepsize () [virtual]

Gets the minimal sampling frequency step size increment value of the device.

Returns

Sampling frequency step size in Hz.

10.46.3.53 GetNumberOfDataBits() virtual uint32_t GetNumberOfDataBits (uint32_t *virtualDevice*, DacqGroupChannelEnumNet *group*, [System::Runtime::InteropServices::Out] int% *numberOfBits*) [virtual]

Get the real number of data bits.

This value may be different from the value returned by GetDataFormat, e.g. in MC_Card the data are shifted 2 bits so the real number is 14 while the data format is 16 bits

10.46.3.54 GetPoti() virtual uint32_t GetPoti (uint32_t *channel*, [System::Runtime::InteropServices::Out] uint32_t% *value*) [virtual]

10.46.3.55 GetResolutionPerDigit() virtual uint32_t GetResolutionPerDigit (uint32_t *virtualDevice*, DacqGroupChannelEnumNet *group*, [System::Runtime::InteropServices::Out] int% *res*, [System::Runtime::InteropServices::Out] int% *resUnit*) [virtual]

10.46.3.56 GetSamplerate() virtual int32_t GetSamplerate (unsigned int *virtualDevice*) [virtual]

Gets the sampling frequency of the device.

Returns

Sampling frequency in Hz.

10.46.3.57 GetVoltageRangeIndex() virtual uint32_t GetVoltageRangeIndex (unsigned int *virtualDevice*) [virtual]

10.46.3.58 GetVoltageRangeInMicroVolt() `virtual int32_t GetVoltageRangeInMicroVolt (unsigned int virtualDevice) [virtual]`

Gets the currently selected voltage range on devices which support multiple voltage ranges.

Returns

The Voltage Range in uV.

10.46.3.59 GetVoltageRangeInMilliVolt() `virtual int32_t GetVoltageRangeInMilliVolt () [virtual]`

Gets the currently selected voltage range on devices which support multiple voltage ranges.

Returns

The rounded Voltage Range in mV.

10.46.3.60 HWInfo() `CHWInfo ^ HWInfo ()`

10.46.3.61 SendStartDacq() [1/2] `virtual void SendStartDacq () [virtual]`

Start sampling.

10.46.3.62 SendStartDacq() [2/2] `virtual void SendStartDacq (int VirtualDacqMap) [virtual]`

Start sampling.

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

10.46.3.63 SendStartStgAndDacq() `virtual void SendStartStgAndDacq (uint32_t trigger_map, int VirtualDacqMap) [virtual]`

Start sampling together with the STG.

Parameters

<i>trigger_map</i>	
<i>VirtualDacqMap</i>	

10.46.3.64 SendStopDacq() [1/2] virtual void SendStopDacq () [virtual]

Stop sampling.

10.46.3.65 SendStopDacq() [2/2] virtual void SendStopDacq (
int *VirtualDacqMap*) [virtual]

Stop sampling.

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

10.46.3.66 SendStopStgAndDacq() virtual void SendStopStgAndDacq (
uint32_t *trigger_map*,
int *VirtualDacqMap*) [virtual]

Stop sampling together with the STG.

Parameters

<i>trigger_map</i>	
--------------------	--

10.46.3.67 SendStopStgAndDacqWithOptions() virtual void SendStopStgAndDacqWithOptions (
uint32_t *trigger_map*,
int *VirtualDacqMap*,
int *options*) [virtual]

Stop sampling together with the STG and options.

Parameters

<i>trigger_map</i>	
--------------------	--

Parameters

<i>options</i>	
----------------	--

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

10.46.3.68 SetDataMode() `virtual void SetDataMode (DataModeEnumNet dataMode, unsigned int virtualDevice) [virtual]`

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

<i>dataMode</i>	DataModeEnumNet enumerates the possible data modes.
<i>virtualDevice</i>	Virtual device to use.

10.46.3.69 SetDigitalSource() `[1/5] void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSource< digitalsourceenum >^ source, int bitnumber_offset)`

Sets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.70 SetDigitalSource() [2/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`DigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.71 SetDigitalSource() [3/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`MEA2100_256DigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.72 SetDigitalSource() [4/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`SCUDigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.73 SetDigitalSource() [5/5] `void SetDigitalSource (`
 DigitalTargetEnumNet *digitaltarget*,
 int32_t *NrChannel*,
 W2100DigitalSourceEnumNet *source*,
 int *bitnumber_offset*)

Sets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

10.46.3.74 SetPoti() `virtual uint32_t SetPoti (`
 uint32_t *channel*,
 uint32_t *value*,
 bool *write_nvram*) [virtual]

10.46.3.75 SetSamplerate() `virtual void SetSamplerate (`
 int32_t *rate*,
 unsigned int *oversample*,
 unsigned int *virtualDevice*) [virtual]

Sets the sampling frequency of the device.

Parameters

<i>rate</i>	Sampling frequency in Hz.
-------------	---------------------------

10.46.3.76 SetSelectedChannels() [1/4] `virtual void SetSelectedChannels (`

```

    array< bool >^ selectedChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    int ChannelsInBlock ) [virtual]

```

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be divided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

10.46.3.77 SetSelectedChannels() [2/4] virtual void SetSelectedChannels (

```

    array< bool >^ selectedChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,

```

```
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

10.46.3.78 SetSelectedChannels() [3/4] virtual void SetSelectedChannels (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock ) [virtual]
```

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be divided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of samples frames the FIFO must acquire before the callback function is called.
------------------	--

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

10.46.3.79 SetSelectedChannels() [4/4] virtual void SetSelectedChannels (

```
int nChannels,
int queuesize,
int threshold,
```

```

SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]

```

10.46.3.80 SetSelectedChannelsQueue() [1/4] virtual void SetSelectedChannelsQueue (
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock) [virtual]

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_↔ ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

10.46.3.81 SetSelectedChannelsQueue() [2/4] virtual void SetSelectedChannelsQueue (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

10.46.3.82 SetSelectedChannelsQueue() [3/4] virtual void SetSelectedChannelsQueue (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use `ChannelBlock_↔ReadFramesDict...` with `handle = 0` to read the data.

When using 32 bit data format, `ChannelsInBlock` is still the number of 16 bit channels per frame, as obtained from `GetChannelsInBlock`, while `nChannels` is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format `nChannels = ChannelsInBlock/2`

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

10.46.3.83 SetSelectedChannelsQueue() [4/4] virtual void SetSelectedChannelsQueue (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

10.46.3.84 SetSelectedData() [1/4] virtual void SetSelectedData (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

10.46.3.85 SetSelectedData() [2/4] virtual void SetSelectedData (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

10.46.3.86 SetSelectedData() [3/4] virtual void SetSelectedData (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

10.46.3.87 SetSelectedData() [4/4] virtual void SetSelectedData (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

10.46.3.88 SetupGroupDacqQueue() void SetupGroupDacqQueue (
 int *queuesize*,
 int *threshold*)

10.46.3.89 SetVoltageRangeByIndex() virtual void SetVoltageRangeByIndex (
 int32_t *voltageRangeIndex*,
 unsigned int *virtualDevice*) [virtual]

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

<i>voltageRangeIndex</i>	Voltage Range to use as index, smaller values are larger voltage ranges.
--------------------------	--

10.46.3.90 SetVoltageRangeInMicroVolt() virtual void SetVoltageRangeInMicroVolt (
 int32_t *voltageRange*,
 unsigned int *virtualDevice*) [virtual]

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

<i>voltageRange</i>	Voltage Range to use in μ V.
---------------------	----------------------------------

This replaces SetVoltageRange, where the value of the range was in mV!

10.46.3.91 StartDacq() [1/4] `virtual void StartDacq () [virtual]`

Start the data acquisition thread and sampling.

10.46.3.92 StartDacq() [2/4] `virtual void StartDacq (
int32_t timeout) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

10.46.3.93 StartDacq() [3/4] `virtual void StartDacq (
int32_t timeout,
int32_t numSubmittedUsbBuffers,
int32_t numUsbBuffers,
int32_t packetsInUrb) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

10.46.3.94 StartDacq() [4/4] `virtual void StartDacq (`
 `int32_t timeout,`
 `int32_t numSubmittedUsbBuffers,`
 `int32_t numUsbBuffers,`
 `int32_t packetsInUrb,`
 `uint32_t virtualDevice) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

10.46.3.95 StartLoop() [1/4] `virtual void StartLoop () [virtual]`

Start the data acquisition thread.

10.46.3.96 StartLoop() [2/4] virtual void StartLoop (
int32_t timeout) [virtual]

Start the data acquisition thread.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

10.46.3.97 StartLoop() [3/4] virtual void StartLoop (
int32_t timeout,
int32_t numSubmittedUsbBuffers,
int32_t numUsbBuffers,
int32_t packetsInUrb) [virtual]

Start the data acquisition thread.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

10.46.3.98 StartLoop() [4/4] virtual void StartLoop (
int32_t timeout,

```

    int32_t numSubmittedUsbBuffers,
    int32_t numUsbBuffers,
    int32_t packetsInUrb,
    uint32_t virtualDevice ) [virtual]

```

Start the data acquisition thread.

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

10.46.3.99 StopDacq() [1/2] virtual void StopDacq () [virtual]

Stop the data acquisition thread and sampling.

10.46.3.100 StopDacq() [2/2] virtual void StopDacq (
 uint32_t virtualDevice) [virtual]

Stop the data acquisition thread and sampling.

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

10.46.3.101 StopLoop() `virtual void StopLoop () [virtual]`

10.46.4 Member Data Documentation

10.46.4.1 Error_Callback_Aquisition_Stopped `const int Error_Callback_Aquisition_Stopped = 0x200 [static]`

10.46.4.2 Error_Callback_Data_lost `const int Error_Callback_Data_lost = 5 [static]`

10.46.4.3 Error_Callback_Frames_Lost `const int Error_Callback_Frames_Lost = 4 [static]`

10.46.4.4 Error_Callback_Packet_Error `const int Error_Callback_Packet_Error = 1 [static]`

10.46.4.5 Error_Callback_Queue_Full `const int Error_Callback_Queue_Full = 0x100 [static]`

10.46.4.6 Error_Callback_RingQueue_Full `const int Error_Callback_RingQueue_Full = 3 [static]`

10.46.5 Property Documentation

10.46.5.1 Samplerate `virtual int Samplerate [get], [set]`


The sampling frequency of the device in Hz.

10.46.6 Event Documentation

10.46.6.1 ChannelDataEvent `OnChannelData^ ChannelDataEvent`10.46.6.2 ErrorEvent `OnError^ ErrorEvent`

10.47 CMcsUsbDeviceStatePushFunctionNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushFunctionNet:

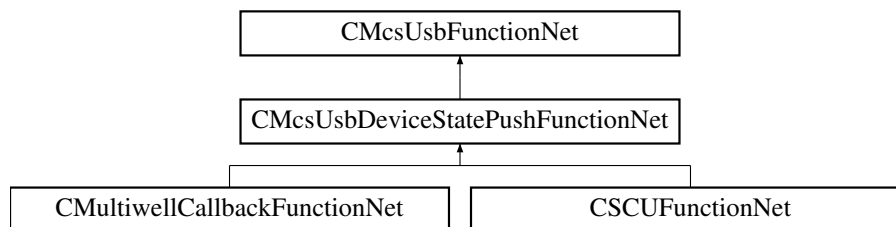

Public Member Functions

- void `TriggerStatus` ()

Protected Member Functions

- `CMcsUsbDeviceStatePushFunctionNet` (`CMcsUsbNet^ mcsusb`, `CMcsUsbFunctionPointerContainer^ pDevice`)

Events

- `OnMcsUsbDeviceState^ McsUsbDeviceStateEvent` [add, remove, raise]

Additional Inherited Members

10.47.1 Constructor & Destructor Documentation

10.47.1.1 CMcsUsbDeviceStatePushFunctionNet()

```

CMcsUsbDeviceStatePushFunctionNet (
    CMcsUsbNet^ mcsusb,
    CMcsUsbFunctionPointerContainer^ pDevice ) [protected]
  
```

10.47.2 Member Function Documentation


10.47.2.1 TriggerStatus() `void TriggerStatus ()`

10.47.3 Event Documentation

10.47.3.1 McsUsbDeviceStateEvent `OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]`

10.48 CMcsUsbDeviceStatePushNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushNet:

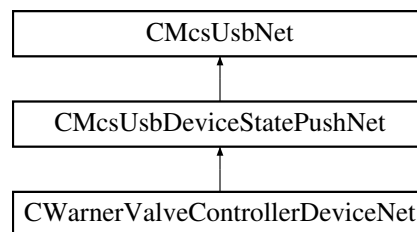

Public Member Functions

- `void TriggerStatus ()`

Protected Member Functions

- `CMcsUsbDeviceStatePushNet (CMcsUsbPointerContainer^ pDevice)`

Events

- `OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add, remove, raise]`

Additional Inherited Members

10.48.1 Constructor & Destructor Documentation

10.48.1.1 CMcsUsbDeviceStatePushNet() `CMcsUsbDeviceStatePushNet (CMcsUsbPointerContainer^ pDevice) [protected]`

10.48.2 Member Function Documentation


10.48.2.1 TriggerStatus() `void TriggerStatus ()`

10.48.3 Event Documentation

10.48.3.1 McsUsbDeviceStateEvent `OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]`

10.49 CMcsUsbFactoryNet Class Reference

Inheritance diagram for CMcsUsbFactoryNet:

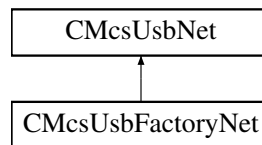

Public Member Functions

- `CMcsUsbFactoryNet ()`
- `~CMcsUsbFactoryNet ()`
- `unsigned int GetNumDestinations ()`
- `String ^ GetDestinationName (unsigned int index)`
- `String ^ GetDestinationName (CFirmwareDestinationNet dest)`
- `void SetDestinationSerialNumber (CFirmwareDestinationNet dest, String^ serialnumber)`
- `String ^ GetDestinationSerialNumber (CFirmwareDestinationNet dest)`
- `CFirmwareDestinationNet GetDestination (unsigned int index)`
- `CFirmwareDestinationNet GetDestination (String^ Key)`
- `unsigned int GetDestinationTargetAddress (CFirmwareDestinationNet destination)`
Gets the target base address for the destination.
- `uint32_t ChangeSerialNumber (String^ serial)`
- `bool LoadUserFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry)`
Sends the DSP Firmware to the MEA21 device.
- `bool LoadUserFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, uint32_t LockMask)`
- `bool UpdateFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange^ deleg, OnUpdateFirmwareProgress^ progress, bool SkipWait)`
Flashes a firmware file to the device.

- bool [UpdateFirmware](#) (String^ FirmwareFile, [CMcsUsbListEntryNet](#)^ listEntry, CFirmwareDestinationNet Dest, [OnUpdateFirmwareStatusChange](#)^ deleg, [OnUpdateFirmwareProgress](#)^ progress, bool SkipWait, unsigned int LockMask)
- bool [UpdateFirmware](#) (String^ FirmwareFile, [CMcsUsbListEntryNet](#)^ listEntry, CFirmwareDestinationNet dest)
- Flashes a firmware file to the device.*
- bool [UpdateFirmware](#) (String^ FirmwareFile, [CMcsUsbListEntryNet](#)^ listEntry, CFirmwareDestinationNet dest, bool SkipWait)
- Flashes a firmware file to the device.*
- bool [UpdateFirmware](#) (String^ FirmwareFile, [CMcsUsbListEntryNet](#)^ listEntry, CFirmwareDestinationNet dest, bool SkipWait, uint32_t LockMask)
- bool [CompareFirmware](#) (String^ FirmwareFile, [CMcsUsbListEntryNet](#)^ listEntry, CFirmwareDestinationNet Dest, [OnUpdateFirmwareStatusChange](#)^ deleg, [OnUpdateFirmwareProgress](#)^ progress, String^ MessagePrefix, unsigned int LockMask, [System::Runtime::InteropServices::Out] String^% ErrorText, [System::Runtime::InteropServices::Out] String^% Protokoll)
- uint32_t [Coldstart](#) (CFirmwareDestinationNet dest)
- int32_t [GetXilinxFlashOffset](#) (CFirmwareDestinationNet dest)
- uint32_t [GetXilinxFlashReadCommand](#) (CFirmwareDestinationNet dest)
- array< uint8_t > ^ [DownloadFirmware](#) (CFirmwareDestinationNet Dest, uint32_t Address, uint32_t length)
- bool [GetUserCodeFromFlash](#) (unsigned int FPGA, unsigned int Address, [System::Runtime::InteropServices::Out] unsigned int% Usercode)
- array< unsigned char > ^ [ReadBlockFromFlash](#) (unsigned int FPGA, unsigned int Address)
- void [ReadBlockFromFlash](#) (unsigned int FPGA, unsigned int Address, array< unsigned char >^ buffer, int position)
- array< unsigned char > ^ [ReadBlockFromIFBGlobalEEPROM](#) (unsigned int Address)
- array< unsigned char > ^ [ReadBlockFromNVMEM](#) (unsigned int FPGA, unsigned int Offset, unsigned int Address)

Static Public Member Functions

- static String ^ [GetDestinationDisplayLabel](#) (String^ RawLabel, CFirmwareDestinationNet dest)
- static String ^ [FindFirmwareVersionMagicInBuffer](#) (array< unsigned char >^ buffer, int length, [System::Runtime::InteropServices::Out] int% position)
- static bool [GetFirmwareVersionFromFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
- Retrives version info from a Firmware update file.*
- static bool [GetFirmwareVersionFromFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version, [System::Runtime::InteropServices::Out] uint32_t% Position)
- static bool [GetFirmwareVersionFromHexFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
- static uint32_t [GetChecksumFromFX3Image](#) (String^ FirmwareFile)
- static uint32_t [GetUSBDeviceIDFromFX3Image](#) (String^ FirmwareFile)
- static bool [GetUserCodeFromBitFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] unsigned int% Usercode)

Static Public Attributes

- static const uint32_t [FX3MCSDDataAddress](#) = 0x40037E00
- static const uint32_t [FX3MCSDDataDeviceIDOffset](#) = 0x4
- static const uint32_t [FX3MCSDDataVersionOffset](#) = 0x8

Additional Inherited Members

10.49.1 Constructor & Destructor Documentation

10.49.1.1 CMcsUsbFactoryNet() `CMcsUsbFactoryNet ()`

10.49.1.2 ~CMcsUsbFactoryNet() `~CMcsUsbFactoryNet ()`

10.49.2 Member Function Documentation

10.49.2.1 ChangeSerialNumber() `uint32_t ChangeSerialNumber (String^ serial)`

10.49.2.2 Coldstart() `uint32_t Coldstart (CFirmwareDestinationNet dest)`

10.49.2.3 CompareFirmware() `bool CompareFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange^ deleg, OnUpdateFirmwareProgress^ progress, String^ MessagePrefix, unsigned int LockMask, [System::Runtime::InteropServices::Out] String^% ErrorText, [System::Runtime::InteropServices::Out] String^% Protokoll)`

10.49.2.4 DownloadFirmware() `array<uint8_t> ^ DownloadFirmware (CFirmwareDestinationNet Dest, uint32_t Address, uint32_t length)`

10.49.2.5 FindFirmwareVersionMagicInBuffer() static String ^ FindFirmwareVersionMagicInBuffer (array< unsigned char >^ *buffer*, int *length*, [System::Runtime::InteropServices::Out] int% *position*) [static]

10.49.2.6 GetChecksumFromFX3Image() static uint32_t GetChecksumFromFX3Image (String^ *FirmwareFile*) [static]

10.49.2.7 GetDestination() [1/2] CFirmwareDestinationNet GetDestination (String^ *Key*)

10.49.2.8 GetDestination() [2/2] CFirmwareDestinationNet GetDestination (unsigned int *index*)

10.49.2.9 GetDestinationDisplayLabel() static String ^ GetDestinationDisplayLabel (String^ *RawLabel*, CFirmwareDestinationNet *dest*) [static]

10.49.2.10 GetDestinationName() [1/2] String ^ GetDestinationName (CFirmwareDestinationNet *dest*)

10.49.2.11 GetDestinationName() [2/2] String ^ GetDestinationName (unsigned int *index*)

10.49.2.12 GetDestinationSerialNumber() String ^ GetDestinationSerialNumber (CFirmwareDestinationNet *dest*)

10.49.2.13 GetDestinationTargetAddress() unsigned int GetDestinationTargetAddress (CFirmwareDestinationNet *destination*)

Gets the target base address for the destination.

Parameters

<i>destination</i>	The destination to be queried.
--------------------	--------------------------------

Returns

The base address as a 32 bit number, only the lower 16 bit represent the address.

10.49.2.14 GetFirmwareVersionFromFile() [1/2] static bool GetFirmwareVersionFromFile (
 String^ *FirmwareFile*,
 [System::Runtime::InteropServices::Out] uint32_t% *Version*) [static]

Retrives version info from a Firmware update file.

10.49.2.15 GetFirmwareVersionFromFile() [2/2] static bool GetFirmwareVersionFromFile (
 String^ *FirmwareFile*,
 [System::Runtime::InteropServices::Out] uint32_t% *Version*,
 [System::Runtime::InteropServices::Out] uint32_t% *Position*) [static]

10.49.2.16 GetFirmwareVersionFromHexFile() static bool GetFirmwareVersionFromHexFile (
 String^ *FirmwareFile*,
 [System::Runtime::InteropServices::Out] uint32_t% *Version*) [static]

10.49.2.17 GetNumDestinations() unsigned int GetNumDestinations ()

10.49.2.18 GetUSBDeviceIDFromFX3Image() static uint32_t GetUSBDeviceIDFromFX3Image (
 String^ *FirmwareFile*) [static]

10.49.2.19 GetUsercodeFromBitFile() static bool GetUsercodeFromBitFile (
 String^ *FirmwareFile*,
 [System::Runtime::InteropServices::Out] unsigned int% *Usercode*) [static]

10.49.2.20 GetUsercodeFromFlash() `bool GetUsercodeFromFlash (unsigned int FPGA, unsigned int Address, [System::Runtime::InteropServices::Out] unsigned int% Usercode)`

10.49.2.21 GetXilinxFlashOffset() `int32_t GetXilinxFlashOffset (CFirmwareDestinationNet dest)`

10.49.2.22 GetXilinxFlashReadCommand() `uint32_t GetXilinxFlashReadCommand (CFirmwareDestinationNet dest)`

10.49.2.23 LoadUserFirmware() [1/2] `bool LoadUserFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry)`

Send the DSP Firmware to the MEA21 device.

Parameters

<i>FirmwareFile</i>	Filename of the DSP Firmware (*.bin) file.
---------------------	--

Parameters

<i>listEntry</i>	Device to use for the connection. See CMcsUsbListNet .
------------------	--

10.49.2.24 LoadUserFirmware() [2/2] `bool LoadUserFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, uint32_t LockMask)`

10.49.2.25 ReadBlockFromFlash() [1/2] `array<unsigned char> ^ ReadBlockFromFlash (unsigned int FPGA, unsigned int Address)`

10.49.2.26 ReadBlockFromFlash() [2/2] void ReadBlockFromFlash (
 unsigned int *FPGA*,
 unsigned int *Address*,
 array< unsigned char >^ *buffer*,
 int *position*)

10.49.2.27 ReadBlockFromIFBGlobalEEProm() array<unsigned char> ^ ReadBlockFromIFBGlobalEEProm (
 unsigned int *Address*)

10.49.2.28 ReadBlockFromNVMEM() array<unsigned char> ^ ReadBlockFromNVMEM (
 unsigned int *FPGA*,
 unsigned int *Offset*,
 unsigned int *Address*)

10.49.2.29 SetDestinationSerialNumber() void SetDestinationSerialNumber (
 CFirmwareDestinationNet *dest*,
 String^ *serialnumber*)

10.49.2.30 UpdateFirmware() [1/5] bool UpdateFirmware (
 String^ *FirmwareFile*,
 CMcsUsbListEntryNet^ *listEntry*,
 CFirmwareDestinationNet *dest*)

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------

Parameters

<i>listEntry</i>	Device to use for the connection.
------------------	-----------------------------------

10.49.2.31 UpdateFirmware() [2/5] bool UpdateFirmware (
 String^ *FirmwareFile*,

```
CMcsUsbListEntryNet^ listEntry,  
CFirmwareDestinationNet dest,  
bool SkipWait )
```

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------

Parameters

<i>listEntry</i>	Device to use for the connection.
------------------	-----------------------------------

10.49.2.32 UpdateFirmware() [3/5] bool UpdateFirmware (
String^ FirmwareFile,
CMcsUsbListEntryNet^ listEntry,
CFirmwareDestinationNet dest,
bool SkipWait,
uint32_t LockMask)

10.49.2.33 UpdateFirmware() [4/5] bool UpdateFirmware (
String^ FirmwareFile,
CMcsUsbListEntryNet^ listEntry,
CFirmwareDestinationNet Dest,
OnUpdateFirmwareStatusChange^ deleg,
OnUpdateFirmwareProgress^ progress,
bool SkipWait)

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------


10.49.2.34 UpdateFirmware() [5/5] bool UpdateFirmware (
String^ FirmwareFile,
CMcsUsbListEntryNet^ listEntry,
CFirmwareDestinationNet Dest,
OnUpdateFirmwareStatusChange^ deleg,

10.49.3 Member Data Documentation

10.49.3.3 FX3MCSDataVersionOffset

```
const uint32_t FX3MCSDataVersionOffset = 0x8 [static]
```

Inheritance diagram for CMcsUsbFunctionNet:

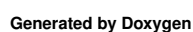

Public Member Functions

- [CMcsUsbFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMcsUsbFunctionNet](#) (void)
- [!CMcsUsbFunctionNet](#) ()
- void [ThrowCUsbExceptionNetOnError](#) (uint32_t status)

Protected Member Functions

- [CMcsUsbFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] mcsusbfunction)

Protected Attributes

- [CMcsUsbNet](#)[^] [m_pMcsUsb](#)
- [CMcsUsbFunction](#) * [m_pMcsUsbFunction](#)

10.50.1 Constructor & Destructor Documentation

10.50.1.1 [CMcsUsbFunctionNet\(\)](#) [1/2] [CMcsUsbFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

10.50.1.2 [~CMcsUsbFunctionNet\(\)](#) virtual [~CMcsUsbFunctionNet](#) (
void) [virtual]

10.50.1.3 [!CMcsUsbFunctionNet\(\)](#) [!CMcsUsbFunctionNet](#) ()

10.50.1.4 [CMcsUsbFunctionNet\(\)](#) [2/2] [CMcsUsbFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *mcsusbfunction*) [protected]

10.50.2 Member Function Documentation

10.50.2.1 [ThrowCUsbExceptionNetOnError\(\)](#) void [ThrowCUsbExceptionNetOnError](#) (
uint32_t *status*)

10.50.3 Member Data Documentation

10.50.3.1 `m_pMcsUsb` [CMcsUsbNet](#) ^ `m_pMcsUsb` [protected]

10.50.3.2 `m_pMcsUsbFunction` `CMcsUsbFunction*` `m_pMcsUsbFunction` [protected]

10.51 CMcsUsbFunctionPointerContainer Class Reference

10.52 CMcsUsbListEntryNet Class Reference

McsUsbListEntryNet identifies a connected device.

Public Member Functions

- `~CMcsUsbListEntryNet` ()
- virtual bool `Equals` (Object^ obj) override
Checks weather two [CMcsUsbListEntryNet](#) represent the same USB device.
- void `SetStringFormat` (String ^ format)
Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.
- virtual String ^ `ToString` () override

Static Public Member Functions

- static `CMcsUsbListEntryNet` ^ `GetEntry` ()
Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.
- static `CMcsUsbListEntryNet` ^ `GetEntry` (DeviceEnumNet McsUsbDevice)
Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.
- static `CMcsUsbListEntryNet` ^ `GetEntry` (DeviceEnumNet McsUsbDevice, unsigned int index)
Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.
- static unsigned int `GetEntryCount` ()
Returns the number of devices connected to the computer.
- static unsigned int `GetEntryCount` (DeviceEnumNet McsUsbDevice)
Returns the number of devices connected to the computer.

Properties

- String^ `Manufacturer` [get]
The Manufacturer ID of the device represented by this [CMcsUsbListEntryNet](#).
- String^ `Product` [get]
The Product ID of the device represented by this [CMcsUsbListEntryNet](#).
- String^ `DeviceName` [get]
The device name of the device represented by this [CMcsUsbListEntryNet](#).
- String^ `SerialNumber` [get]
The serial number of the device represented by this [CMcsUsbListEntryNet](#).
- String^ `HwVersion` [get]
The hardware revision of the device represented by this [CMcsUsbListEntryNet](#).
- `DeviceldNet`^ `Deviceld` [get]

10.52.1 Detailed Description

CMcsUsbListEntryNet identifies a connected device.

10.52.2 Constructor & Destructor Documentation

10.52.2.1 `~CMcsUsbListEntryNet()` `~CMcsUsbListEntryNet ()`

10.52.3 Member Function Documentation

10.52.3.1 `Equals()` `virtual bool Equals (`
`Object^ obj) [override], [virtual]`

Checks weather two [CMcsUsbListEntryNet](#) represent the same USB device.

Parameters

<i>obj</i>	The CMcsUsbListEntryNet to compare with.
------------	--

10.52.3.2 `GetEntry()` `[1/3] static CMcsUsbListEntryNet ^ GetEntry () [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

10.52.3.3 `GetEntry()` `[2/3] static CMcsUsbListEntryNet ^ GetEntry (`
`DeviceEnumNet McsUsbDevice) [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

10.52.3.4 GetEntry() [3/3] `static CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice, unsigned int index) [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Parameters

<i>index</i>	number of the entry to use.
--------------	-----------------------------

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

10.52.3.5 GetEntryCount() [1/2] `static unsigned int GetEntryCount () [static]`

Returns the number of devices connected to the computer.

Returns

The number of devices.

10.52.3.6 GetEntryCount() [2/2] `static unsigned int GetEntryCount (DeviceEnumNet McsUsbDevice) [static]`

Returns the number of devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Returns

The number of devices.

10.52.3.7 SetStringFormat() `void SetStringFormat (String ^ format)`

Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

<i>format</i>	A String containing the format template. Default is "%N (%S)".
---------------	--

10.52.3.8 ToString() `virtual String ^ ToString () [override], [virtual]`

10.52.4 Property Documentation

10.52.4.1 DeviceId `DeviceIdNet^ DeviceId [get]`

10.52.4.2 DeviceName `String^ DeviceName [get]`

The device name of the device represented by this [CMcsUsbListEntryNet](#).

10.52.4.3 HwVersion `String^ HwVersion [get]`

The hardware revision of the device represented by this [CMcsUsbListEntryNet](#).

10.52.4.4 Manufacturer `String^ Manufacturer [get]`

The Manufacturer ID of the device represented by this [CMcsUsbListEntryNet](#).

10.52.4.5 Product `String^ Product [get]`

The Product ID of the device represented by this [CMcsUsbListEntryNet](#).

10.52.4.6 SerialNumber `String^ SerialNumber [get]`

The serial number of the device represented by this [CMcsUsbListEntryNet](#).

10.53 CMcsUsbListNet Class Reference

Class to handle a list of connected MCS USB devices.

Public Member Functions

- [CMcsUsbListNet](#) (DeviceEnumNet McsUsbDevice)
Initializes a new instance of [CMcsUsbListNet](#) class.
- [CMcsUsbListNet](#) (array< [DeviceIdNet](#) ^> ^ DeviceIdList)
Initializes a new instance of [CMcsUsbListNet](#) class.
- [~CMcsUsbListNet](#) ()
Destructor: called by [Dispose\(\)](#)
- [!CMcsUsbListNet](#) ()
Finalizer: called by GC before collecting
- void [SetStringFormat](#) (String ^ format)
Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.
- uint32_t [GetNumberOfDevices](#) ()
Gets the number of devices currently in the list.
- [CMcsUsbListEntryNet](#) ^ [GetUsbListEntry](#) (unsigned int index)
Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.
- array< [CMcsUsbListEntryNet](#) ^> ^ [GetUsbListEntries](#) ()
Returns all entries from the list of USB Devices connected to the computer.
- bool [IsDeviceTypeOf](#) ([CMcsUsbListEntryNet](#) ^ entry, DeviceEnumNet McsUsbDevice)

Properties

- uint32_t [Count](#) [get]
Gets the number of devices currently in the list.

Events

- [OnDeviceArrivalRemoval](#)^ [DeviceArrival](#)
- [OnDeviceArrivalRemoval](#)^ [DeviceRemoval](#)

10.53.1 Detailed Description

Class to handle a list of connected MCS USB devices.

10.53.2 Constructor & Destructor Documentation

10.53.2.1 CMcsUsbListNet() [1/2] [CMcsUsbListNet](#) (
 [DeviceEnumNet](#) *McsUsbDevice*)

Initializes a new instance of [CMcsUsbListNet](#) class.

10.53.2.2 CMcsUsbListNet() [2/2] [CMcsUsbListNet](#) (
 [array](#)< [DeviceIdNet](#)^>^ [DeviceIdList](#))

Initializes a new instance of [CMcsUsbListNet](#) class.

10.53.2.3 ~CMcsUsbListNet() [~CMcsUsbListNet](#) ()

Destructor: called by Dispose()

10.53.2.4 ~!CMcsUsbListNet() [!CMcsUsbListNet](#) ()

Finalizer: called by GC before collecting

10.53.3 Member Function Documentation

10.53.3.1 GetNumberOfDevices() [uint32_t](#) [GetNumberOfDevices](#) ()

Gets the number of devices currently in the list.

Returns

The number of devices currently in the list.

10.53.3.2 GetUsbListEntries() [array](#)<[CMcsUsbListEntryNet](#)^> ^ [GetUsbListEntries](#) ()

Returns all entries from the list of USB Devices connected to the computer.

10.53.3.3 GetUsbListEntry() [CMcsUsbListEntryNet](#) ^ [GetUsbListEntry](#) (
 [unsigned int](#) *index*)

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Parameters

<i>index</i>	number of the entry to use.
--------------	-----------------------------

10.53.3.4 IsDeviceTypeOf() `bool IsDeviceTypeOf (`
 `CMcsUsbListEntryNet^ entry,`
 `DeviceEnumNet McsUsbDevice)`

10.53.3.5 SetStringFormat() `void SetStringFormat (`
 `String ^ format)`

Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

<i>format</i>	A String containing the format template. Default is "%N (%S)".
---------------	--

10.53.4 Property Documentation

10.53.4.1 Count `uint32_t Count [get]`

Gets the number of devices currently in the list.

10.53.5 Event Documentation


10.53.5.1 DeviceArrival `OnDeviceArrivalRemoval^ DeviceArrival`

10.53.5.2 DeviceRemoval `OnDeviceArrivalRemoval^ DeviceRemoval`

10.54 CMcsUsbNet Class Reference

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

Inheritance diagram for CMcsUsbNet:

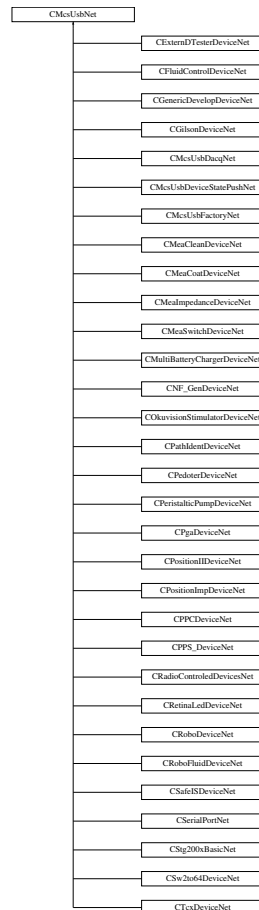

Public Member Functions

- [CMcsUsbNet](#) ()
Initializes a new instance of the base class to handle MCS USB devices.
- [CMcsUsbNet](#) (McsBusTypeEnumNet bustype)
Initializes a new instance of the base class to handle MCS USB devices.
- virtual [~CMcsUsbNet](#) ()
- [!CMcsUsbNet](#) ()
- DeviceEnumNet [GetDeviceEnum](#) ()
- virtual uint32_t [Connect](#) ()
Opens a connection to the device.
- virtual uint32_t [Connect](#) (unsigned int LockMask)
Opens a connection to the device.
- virtual uint32_t [Connect](#) ([CMcsUsbListEntryNet](#)[^] entry)
Opens a connection to the device.
- virtual uint32_t [Connect](#) ([CMcsUsbListEntryNet](#)[^] entry, unsigned int LockMask)
Opens a connection to the device.

- virtual uint32_t [GetStatus](#) ([System::Runtime::InteropServices::Out]uint32_t% iStatus)
- virtual bool [IsConnected](#) ()
Check if a device is Connected.
- virtual void [Disconnect](#) ()
Disconnect from a device.
- [CMcsUsbListEntryNet](#) ^ [GetUsbListEntry](#) ()
- virtual String ^ [GetSerialNumber](#) ()
Query the Serial Number of the device.
- [DriverVersionNet](#) ^ [GetVersion](#) ()
- [DriverVersionNet](#) ^ [GetVersion](#) (CFirmwareDestinationNet dest)
- [DeviceldNet](#) ^ [GetDeviceld](#) ()
- uint32_t [GetIdent](#) ([System::Runtime::InteropServices::Out]String^% Answer)
- void [MultibootSelectImage](#) (unsigned int sector)
Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.
- String ^ [MultibootGetImageld](#) (unsigned int sector)
Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.
- uint32_t [MultibootGetCypressImageld](#) (unsigned int sector)
Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.
- uint32_t [MultibootGetSelectedImage](#) ()
Gets sector index of selected FPGA boot image on IFB
- uint32_t [GetMea21UsbPort](#) ()
Gets the USB port if an IFB that is used by this connection
- HeadstageIdEnumNet [GetHeadstageID](#) (uint32_t headstage)
Gets the ID of a connected headstage.
- bool [GetHeadstagePresent](#) (uint32_t headstage)
queries whether a headstage is present
- bool [GetHeadstageActive](#) (uint32_t headstage)
queries whether a headstage is active
- void [RescanHeadstage](#) (uint32_t headstage)
rescans and activates a headstage
- array< BYTE > ^ [GetSoftwareKey](#) (unsigned int index)
- void [SetSoftwareKey](#) (unsigned int index, array< BYTE >^ buffer)
- void [RemoveSoftwareKey](#) (unsigned int index)
- void [AddSoftwareKey](#) (String^ key)
- bool [EmptyKey](#) (String^ key)
- bool [ValidKey](#) (String^ key, [System::Runtime::InteropServices::Out]String^% serial_number)
- bool [ValidKey](#) (String^ key, uint8_t ProgrammID, uint8_t majorversion, [System::Runtime::InteropServices::Out]String^% serial_number)
- bool [HasSoftwareKey](#) (uint8_t ProgrammID, uint8_t majorversion)
- bool [HasSoftwareKey](#) (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- String ^ [GetSoftwareKeyString](#) (uint8_t ProgrammID, uint8_t majorversion)
- String ^ [GetSoftwareKeyString](#) (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- bool [IsDeviceHighSpeedCapable](#) ()
- bool [IsDeviceHighSpeed](#) ()
- McsUsbSpeedEnumNet [GetDeviceCapableSpeed](#) ()
- McsUsbSpeedEnumNet [GetDeviceSpeed](#) ()
Query the Connection Speed of the device.
- unsigned int [TxnTestMemoryWrite](#) (unsigned short index)
- unsigned int [TxnTestMemoryReadAndCheck](#) (unsigned short index)
- void [TxnSetSerialNumber](#) (unsigned int number)
- unsigned int [TxnGetSerialNumber](#) ()
- unsigned int [ReadRegister](#) (unsigned int reg)

- array< uint32_t > ^ [ReadRegister](#) (unsigned int reg, int length)
- unsigned int [ReadRegister32](#) (unsigned int adr)
- unsigned int [ReadRegisterTimeSlot](#) (unsigned int reg, int TimeSlot)
- void [WriteRegister](#) (unsigned int reg, unsigned int value)
- void [WriteRegisterValue](#) (unsigned int reg, unsigned int value)
- void [WriteRegister32](#) (unsigned int adr, unsigned int value)
- void [WriteRegister](#) (unsigned int reg, array< unsigned int >^ values)
- void [WriteRegisterArray](#) (unsigned int reg, array< unsigned int >^ values)
- void [WriteRegisterTimeSlot](#) (unsigned int reg, unsigned int value, int TimeSlot)
- void [WriteRegisterTimeSlot](#) (unsigned int reg, array< unsigned int >^ values, int TimeSlot)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t% DMA_value)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- unsigned int [GetLastUSBError](#) ()
- void [ThrowCUsbExceptionNetOnError](#) (uint32_t status)
- uint32_t [GetDeviceCannotStallOutRequests](#) ()
- String ^ [GetHardwareRevision](#) ()
- unsigned int [GetFirmwareVersion](#) (CFirmwareDestinationNet destination)
Gets the firmware version for the destination.
- uint8_t [GetNumConfigurations](#) ()
- uint8_t [GetConfiguration](#) ()
- void [SetConfiguration](#) (uint8_t config)
- uint32_t [GetDeviceRootHubVendorID](#) ()
Gets the Vendor ID of the USB root hub the device is connected to.
- UsbVendorIdEnumNet [GetDeviceRootHubVendorEnum](#) ()
Gets the Vendor ID of the USB root hub the device is connected to.
- String ^ [GetDeviceRootHubVendorName](#) ()
Gets the Vendor Name of the USB root hub the device is connected to.
- void [EnableExceptions](#) (bool enable)
Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.
- bool [IsExceptionsEnabled](#) ()
- uint32_t [GetStatusOfLastCommand](#) ()
Gets the status of the last call to the McsUsb Library.
- void [AssociateToThis](#) (CMcsUsbNet^ device)

Static Public Member Functions

- static String ^ [GetErrorText](#) (unsigned int Status)
Gets the error text string that belongs to a status number.

Static Public Attributes

- static const uint32_t [Status_Crc](#) = (0xE0100001L)
- static const uint32_t [Status_Btstuff](#) = (0xE0100002L)
- static const uint32_t [Status_DataToggleMismatch](#) = (0xE0100003L)
- static const uint32_t [Status_Stall](#) = (0xE0100004L)
- static const uint32_t [Status_DevNotResponding](#) = (0xE0100005L)
- static const uint32_t [Status_PidCheckFailure](#) = (0xE0100006L)
- static const uint32_t [Status_UnexpectedPid](#) = (0xE0100007L)
- static const uint32_t [Status_DataOverrun](#) = (0xE0100008L)
- static const uint32_t [Status_DataUnderrun](#) = (0xE0100009L)
- static const uint32_t [Status_BufferOverrun](#) = (0xE010000CL)
- static const uint32_t [Status_BufferUnderrun](#) = (0xE010000DL)
- static const uint32_t [Status_NotAccessed](#) = (0xE010000FL)
- static const uint32_t [Status_Fifo](#) = (0xE0100010L)
- static const uint32_t [Status_EndpointHalted](#) = (0xE0100030L)
- static const uint32_t [Status_NoMemory](#) = (0xE0100100L)
- static const uint32_t [Status_InvalidUrbFunction](#) = (0xE0100200L)
- static const uint32_t [Status_InvalidParameter](#) = (0xE0100300L)
- static const uint32_t [Status_ErrorBusy](#) = (0xE0100400L)
- static const uint32_t [Status_RequestFailed](#) = (0xE0100500L)
- static const uint32_t [Status_InvalidPipeHandle](#) = (0xE0100600L)
- static const uint32_t [Status_NoBandwidth](#) = (0xE0100700L)
- static const uint32_t [Status_InternalHcError](#) = (0xE0100800L)
- static const uint32_t [Status_ErrorShortTransfer](#) = (0xE0100900L)
- static const uint32_t [Status_BadStartFrame](#) = (0xE0100A00L)
- static const uint32_t [Status_IsochRequestFailed](#) = (0xE0100B00L)
- static const uint32_t [Status_FrameControlOwned](#) = (0xE0100C00L)
- static const uint32_t [Status_ControlNotOwned](#) = (0xE0100D00L)
- static const uint32_t [Status_Canceled](#) = (0xE0110000L)
- static const uint32_t [Status_Canceling](#) = (0xE0120000L)
- static const uint32_t [Status_AlreadyConfigured](#) = (0xE0110001L)
- static const uint32_t [Status_Unconfigured](#) = (0xE0110002L)
- static const uint32_t [Status_NoSuchDevice](#) = (0xE01F0002L)
- static const uint32_t [Status_DeviceNotFound](#) = (0xE01F0003L)
- static const uint32_t [Status_NotSupported](#) = (0xE01F0005L)
- static const uint32_t [Status_IoPending](#) = (0xE01F0006L)
- static const uint32_t [Status_IoTimeout](#) = (0xE01F0007L)
- static const uint32_t [Status_DeviceRemoved](#) = (0xE01F0008L)
- static const uint32_t [Status_PipeNotLinked](#) = (0xE01F0009L)
- static const uint32_t [Status_ConnectedPipes](#) = (0xE01F000AL)
- static const uint32_t [Status_DeviceLocked](#) = (0xE01F0010L)
- static const uint32_t [Status_RequestMutexTimeout](#) = (0xE01F0020L)
- static const uint32_t [Status_RequestMutexFailed](#) = (0xE01F0021L)
- static const uint32_t [Status_LastUsbErrorMismatch](#) = (0xE01F0022L)
- static const uint32_t [WPAError_ScanningIsPending](#) = ((0xA0220000L) | 0x0036)

Properties

- virtual String[^] [SerialNumber](#) [get]

10.54.1 Detailed Description

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

10.54.2 Constructor & Destructor Documentation

10.54.2.1 CMcsUsbNet() [1/2] `CMcsUsbNet ()`

Initializes a new instance of the base class to handle MCS USB devices.

10.54.2.2 CMcsUsbNet() [2/2] `CMcsUsbNet (McsBusTypeEnumNet bustype)`

Initializes a new instance of the base class to handle MCS USB devices.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

10.54.2.3 ~CMcsUsbNet() `virtual ~CMcsUsbNet () [virtual]`

10.54.2.4 "!CMcsUsbNet() `!CMcsUsbNet ()`

10.54.3 Member Function Documentation

10.54.3.1 AddSoftwareKey() `void AddSoftwareKey (String^ key)`

10.54.3.2 AssociateToThis() `void AssociateToThis (CMcsUsbNet^ device)`

10.54.3.3 Connect() [1/4] `virtual uint32_t Connect () [virtual]`

Opens a connection to the device.

Returns

Error Status. 0 on success.

10.54.3.4 Connect() [2/4] `virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry) [virtual]`

Opens a connection to the device.

Parameters

<i>entry</i>	The Device List Entry for the device to be connected.
--------------	---

Returns

Error Status. 0 on success.

10.54.3.5 Connect() [3/4] `virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry,
unsigned int LockMask) [virtual]`

Opens a connection to the device.

Parameters

<i>entry</i>	The Device List Entry for the device to be connected.
<i>LockMask</i>	The Lock Mask for this connection.

Returns

Error Status. 0 on success.

10.54.3.6 Connect() [4/4] `virtual uint32_t Connect (
unsigned int LockMask) [virtual]`

Opens a connection to the device.

Parameters

<i>LockMask</i>	The Lock Mask for this connection.
-----------------	------------------------------------

Returns

Error Status. 0 on success.

10.54.3.7 Disconnect() `virtual void Disconnect () [virtual]`

Disconnect from a device.

10.54.3.8 EmptyKey() `bool EmptyKey (String^ key)`

10.54.3.9 EnableExceptions() `void EnableExceptions (bool enable)`

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the `GetStatusOfLastCommand` call instead.

Parameters

<i>enable</i>	True to enable Exceptions, False to disable.
---------------	--

10.54.3.10 EraseEepromRegisterPreconfig() [1/3] `void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg)`

10.54.3.11 EraseEepromRegisterPreconfig() [2/3] `void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize)`

10.54.3.12 EraseEepromRegisterPreconfig() [3/3] `void EraseEepromRegisterPreconfig (`
 `uint32_t EEPROMBase,`
 `uint32_t DMA_reg,`
 `uint32_t EEPROMSize,`
 `uint32_t EepromStartAddress)`

10.54.3.13 GetConfiguration() `uint8_t GetConfiguration ()`

10.54.3.14 GetDeviceCannotStallOutRequests() `uint32_t GetDeviceCannotStallOutRequests ()`

10.54.3.15 GetDeviceCapableSpeed() `McsUsbSpeedEnumNet GetDeviceCapableSpeed ()`

10.54.3.16 GetDeviceEnum() `DeviceEnumNet GetDeviceEnum ()`

10.54.3.17 GetDeviceId() `DeviceIdNet ^ GetDeviceId ()`

10.54.3.18 GetDeviceRootHubVendorEnum() `UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ()`

Gets the Vendor ID of the USB root hub the device is connected to.

Returns

An enum which enumerates the PCI Vendor ID.

10.54.3.19 GetDeviceRootHubVendorID() `uint32_t GetDeviceRootHubVendorID ()`

Gets the Vendor ID of the USB root hub the device is connected to.

Returns

The PCI Vendor ID, 0x8086 for Intel, 0x1912 for Renesas, 0x1b21 for ASMedia.

10.54.3.20 GetDeviceRootHubVendorName() `String ^ GetDeviceRootHubVendorName ()`

Gets the Vendor Name of the USB root hub the device is connected to.

Returns

The PCI Vendor Name, either "Intel", "Renesas", "ASMedia" or "unknown".

10.54.3.21 GetDeviceSpeed() `McsUsbSpeedEnumNet GetDeviceSpeed ()`

Query the Connection Speed of the device.

Returns

0 for Low-Speed, 1 for Full-Speed, 2 for High-Speed and 3 for SuperSpeed.

10.54.3.22 GetErrorText() `static String ^ GetErrorText (unsigned int Status) [static]`

Gets the error text string that belongs to a status number.

Parameters

<i>Status</i>	The status number you want the text for.
---------------	--

Returns

The error text string that belongs to the status number.

10.54.3.23 GetFirmwareVersion() `unsigned int GetFirmwareVersion (CFirmwareDestinationNet destination)`

Gets the firmware version for the destination.

Parameters

<i>destination</i>	The destination to be queried.
--------------------	--------------------------------

Returns

The firmware version as a 32 bit number, the upper 16 bit contain the major version number, the lower 16 bit the minor version number.

10.54.3.24 GetHardwareRevision() `String ^ GetHardwareRevision ()`

10.54.3.25 GetHeadstageActive() `bool GetHeadstageActive (
 uint32_t headstage)`

queries whether a headstage is active

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

true if the headstage is active

10.54.3.26 GetHeadstageID() `HeadstageIdEnumNet GetHeadstageID (
 uint32_t headstage)`

Gets the ID of a connected headstage.

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

enumerated Headstage ID

10.54.3.27 GetHeadstagePresent() `bool GetHeadstagePresent (
 uint32_t headstage)`

queries whether a headstage is present

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

true if the headstage is present

10.54.3.28 GetIdent() `uint32_t GetIdent ([System::Runtime::InteropServices::Out] String^% Answer)`

10.54.3.29 GetLastUSBError() `unsigned int GetLastUSBError ()`

10.54.3.30 GetMea21UsbPort() `uint32_t GetMea21UsbPort ()`

Gets the USB port if an IFB that is used by this connection

Returns

number of used port; range: 0..1

10.54.3.31 GetNumConfigurations() `uint8_t GetNumConfigurations ()`

10.54.3.32 GetSerialNumber() `virtual String ^ GetSerialNumber () [virtual]`

Query the Serial Number of the device.

Returns

The Serial Number.

10.54.3.33 GetSoftwareKey() `array<BYTE> ^ GetSoftwareKey (unsigned int index)`

10.54.3.34 GetSoftwareKeyString() [1/2] `String ^ GetSoftwareKeyString (SoftwareKeyProgramIdsNet::ProgramIdsNet ProgrammID, uint8_t majorversion)`

10.54.3.35 GetSoftwareKeyString() [2/2] `String ^ GetSoftwareKeyString (uint8_t ProgrammID, uint8_t majorversion)`

10.54.3.36 GetStatus() virtual uint32_t GetStatus (
 [System::Runtime::InteropServices::Out] uint32_t% iStatus) [virtual]

10.54.3.37 GetStatusOfLastCommand() uint32_t GetStatusOfLastCommand ()

Gets the status of the last call to the McsUsb Library.

Returns

The Error Status of the last McsUsb command. 0 on success.

10.54.3.38 GetUsbListEntry() CMcsUsbListEntryNet ^ GetUsbListEntry ()

10.54.3.39 GetVersion() [1/2] DriverVersionNet ^ GetVersion ()

10.54.3.40 GetVersion() [2/2] DriverVersionNet ^ GetVersion (
 CFirmwareDestinationNet dest)

10.54.3.41 HasSoftwareKey() [1/2] bool HasSoftwareKey (
 SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
 uint8_t majorversion)

10.54.3.42 HasSoftwareKey() [2/2] bool HasSoftwareKey (
 uint8_t ProgrammID,
 uint8_t majorversion)

10.54.3.43 IsConnected() virtual bool IsConnected () [virtual]

Check if a device is Connected.

Returns

true if the device is connected.

10.54.3.44 IsDeviceHighSpeed() `bool IsDeviceHighSpeed ()`

10.54.3.45 IsDeviceHighSpeedCapable() `bool IsDeviceHighSpeedCapable ()`

10.54.3.46 IsExceptionsEnabled() `bool IsExceptionsEnabled ()`

10.54.3.47 MultibootGetCypressImageId() `uint32_t MultibootGetCypressImageId (unsigned int sector)`

Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.

Returns

The magic ident code of the image.

10.54.3.48 MultibootGetImageId() `String ^ MultibootGetImageId (unsigned int sector)`

Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.

Returns

The magic ident code of the image.

10.54.3.49 MultibootGetSelectedImage() `uint32_t MultibootGetSelectedImage ()`

Gets sector index of selected FPGA boot image on IFB

Returns

Sector index of image; range: 0..2

10.54.3.50 MultibootSelectImage() `void MultibootSelectImage (unsigned int sector)`

Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.

Returns

Throws exception on error.

10.54.3.51 ReadEepromRegisterPreconfig() [1/3] `bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out] uint32_t% DMA_value)`

10.54.3.52 ReadEepromRegisterPreconfig() [2/3] `bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out] uint32_t% DMA_value, uint32_t EEPROMSize)`

10.54.3.53 ReadEepromRegisterPreconfig() [3/3] `bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out] uint32_t% DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)`

10.54.3.54 ReadRegister() [1/2] `unsigned int ReadRegister (unsigned int reg)`

10.54.3.55 ReadRegister() [2/2] `array<uint32_t> ^ ReadRegister (unsigned int reg, int length)`

10.54.3.56 ReadRegister32() `unsigned int ReadRegister32 (unsigned int adr)`

10.54.3.57 ReadRegisterTimeSlot() unsigned int ReadRegisterTimeSlot (
 unsigned int *reg*,
 int *TimeSlot*)

10.54.3.58 RemoveSoftwareKey() void RemoveSoftwareKey (
 unsigned int *index*)

10.54.3.59 RescanHeadstage() void RescanHeadstage (
 uint32_t *headstage*)

rescans and activates a headstage

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

10.54.3.60 SetConfiguration() void SetConfiguration (
 uint8_t *config*)

10.54.3.61 SetSoftwareKey() void SetSoftwareKey (
 unsigned int *index*,
 array< BYTE >^ *buffer*)

10.54.3.62 ThrowCUsbExceptionNetOnError() void ThrowCUsbExceptionNetOnError (
 uint32_t *status*)

10.54.3.63 TxnGetSerialNumber() unsigned int TxnGetSerialNumber ()

10.54.3.64 TxnSetSerialNumber() void TxnSetSerialNumber (
 unsigned int *number*)

10.54.3.65 TxnTestMemoryReadAndCheck() unsigned int TxnTestMemoryReadAndCheck (
 unsigned short *index*)

10.54.3.66 TxnTestMemoryWrite() unsigned int TxnTestMemoryWrite (
 unsigned short *index*)

10.54.3.67 ValidKey() [1/2] bool ValidKey (
 String^ *key*,
 [System::Runtime::InteropServices::Out] String^% *serial_number*)

10.54.3.68 ValidKey() [2/2] bool ValidKey (
 String^ *key*,
 uint8_t *ProgrammID*,
 uint8_t *majorversion*,
 [System::Runtime::InteropServices::Out] String^% *serial_number*)

10.54.3.69 WriteEepromRegisterPreconfig() [1/3] void WriteEepromRegisterPreconfig (
 uint32_t *EEPROMBase*,
 uint32_t *DMA_reg*,
 uint32_t *DMA_value*)

10.54.3.70 WriteEepromRegisterPreconfig() [2/3] void WriteEepromRegisterPreconfig (
 uint32_t *EEPROMBase*,
 uint32_t *DMA_reg*,
 uint32_t *DMA_value*,
 uint32_t *EEPROMSize*)

10.54.3.71 WriteEepromRegisterPreconfig() [3/3] void WriteEepromRegisterPreconfig (
 uint32_t *EEPROMBase*,
 uint32_t *DMA_reg*,
 uint32_t *DMA_value*,
 uint32_t *EEPROMSize*,
 uint32_t *EepromStartAddress*)

10.54.3.72 WriteRegister() [1/2] void WriteRegister (
 unsigned int *reg*,
 array< unsigned int >^ *values*)

10.54.3.73 WriteRegister() [2/2] void WriteRegister (
 unsigned int *reg*,
 unsigned int *value*)

10.54.3.74 WriteRegister32() void WriteRegister32 (
 unsigned int *adr*,
 unsigned int *value*)

10.54.3.75 WriteRegisterArray() void WriteRegisterArray (
 unsigned int *reg*,
 array< unsigned int >^ *values*)

10.54.3.76 WriteRegisterTimeSlot() [1/2] void WriteRegisterTimeSlot (
 unsigned int *reg*,
 array< unsigned int >^ *values*,
 int *TimeSlot*)

10.54.3.77 WriteRegisterTimeSlot() [2/2] void WriteRegisterTimeSlot (
 unsigned int *reg*,
 unsigned int *value*,
 int *TimeSlot*)

10.54.3.78 WriteRegisterValue() void WriteRegisterValue (
 unsigned int *reg*,
 unsigned int *value*)

10.54.4 Member Data Documentation

10.54.4.1 Status_AlreadyConfigured `const uint32_t Status_AlreadyConfigured = (0xE0110001L)`
[static]

10.54.4.2 Status_BadStartFrame `const uint32_t Status_BadStartFrame = (0xE0100A00L)` [static]

10.54.4.3 Status_Btstuff `const uint32_t Status_Btstuff = (0xE0100002L)` [static]

10.54.4.4 Status_BufferOverrun `const uint32_t Status_BufferOverrun = (0xE010000CL)` [static]

10.54.4.5 Status_BufferUnderrun `const uint32_t Status_BufferUnderrun = (0xE010000DL)` [static]

10.54.4.6 Status_Canceled `const uint32_t Status_Canceled = (0xE0110000L)` [static]

10.54.4.7 Status_Canceling `const uint32_t Status_Canceling = (0xE0120000L)` [static]

10.54.4.8 Status_ConnectedPipes `const uint32_t Status_ConnectedPipes = (0xE01F000AL)` [static]

10.54.4.9 Status_ControlNotOwned `const uint32_t Status_ControlNotOwned = (0xE0100D00L)` [static]

10.54.4.10 Status_Crc `const uint32_t Status_Crc = (0xE0100001L)` [static]

10.54.4.11 Status_DataOverrun `const uint32_t Status_DataOverrun = (0xE0100008L)` [static]

10.54.4.12 Status_DataToggleMismatch `const uint32_t Status_DataToggleMismatch = (0xE0100003L)`
[static]

10.54.4.13 Status_DataUnderrun `const uint32_t Status_DataUnderrun = (0xE0100009L)` [static]

10.54.4.14 Status_DeviceLocked `const uint32_t Status_DeviceLocked = (0xE01F0010L)` [static]

10.54.4.15 Status_DeviceNotFound `const uint32_t Status_DeviceNotFound = (0xE01F0003L)` [static]

10.54.4.16 Status_DeviceRemoved `const uint32_t Status_DeviceRemoved = (0xE01F0008L)` [static]

10.54.4.17 Status_DevNotResponding `const uint32_t Status_DevNotResponding = (0xE0100005L)`
[static]

10.54.4.18 Status_EndpointHalted `const uint32_t Status_EndpointHalted = (0xE0100030L)` [static]

10.54.4.19 Status_ErrorBusy `const uint32_t Status_ErrorBusy = (0xE0100400L)` [static]

10.54.4.20 Status_ErrorShortTransfer `const uint32_t Status_ErrorShortTransfer = (0xE0100900L)`
[static]

10.54.4.21 Status_Fifo `const uint32_t Status_Fifo = (0xE0100010L)` [static]

10.54.4.22 Status_FrameControlOwned `const uint32_t Status_FrameControlOwned = (0xE0100C00L)`
[static]

10.54.4.23 Status_InternalHcError `const uint32_t Status_InternalHcError = (0xE0100800L) [static]`

10.54.4.24 Status_InvalidParameter `const uint32_t Status_InvalidParameter = (0xE0100300L) [static]`

10.54.4.25 Status_InvalidPipeHandle `const uint32_t Status_InvalidPipeHandle = (0xE0100600L) [static]`

10.54.4.26 Status_InvalidUrbFunction `const uint32_t Status_InvalidUrbFunction = (0xE0100200L) [static]`

10.54.4.27 Status_IoPending `const uint32_t Status_IoPending = (0xE01F0006L) [static]`

10.54.4.28 Status_IoTimeout `const uint32_t Status_IoTimeout = (0xE01F0007L) [static]`

10.54.4.29 Status_IsochRequestFailed `const uint32_t Status_IsochRequestFailed = (0xE0100B00L) [static]`

10.54.4.30 Status_LastUsbErrorMismatch `const uint32_t Status_LastUsbErrorMismatch = (0xE01F0022L) [static]`

10.54.4.31 Status_NoBandwidth `const uint32_t Status_NoBandwidth = (0xE0100700L) [static]`

10.54.4.32 Status_NoMemory `const uint32_t Status_NoMemory = (0xE0100100L) [static]`

10.54.4.33 Status_NoSuchDevice `const uint32_t Status_NoSuchDevice = (0xE01F0002L) [static]`

10.54.4.34 Status_NotAccessed `const uint32_t Status_NotAccessed = (0xE010000FL) [static]`

10.54.4.35 Status_NotSupported `const uint32_t Status_NotSupported = (0xE01F0005L) [static]`

10.54.4.36 Status_PidCheckFailure `const uint32_t Status_PidCheckFailure = (0xE0100006L) [static]`

10.54.4.37 Status_PipeNotLinked `const uint32_t Status_PipeNotLinked = (0xE01F0009L) [static]`

10.54.4.38 Status_RequestFailed `const uint32_t Status_RequestFailed = (0xE0100500L) [static]`

10.54.4.39 Status_RequestMutexFailed `const uint32_t Status_RequestMutexFailed = (0xE01F0021L) [static]`

10.54.4.40 Status_RequestMutexTimeout `const uint32_t Status_RequestMutexTimeout = (0xE01F0020L) [static]`

10.54.4.41 Status_Stall `const uint32_t Status_Stall = (0xE0100004L) [static]`

10.54.4.42 Status_Unconfigured `const uint32_t Status_Unconfigured = (0xE0110002L) [static]`

10.54.4.43 Status_UnexpectedPid `const uint32_t Status_UnexpectedPid = (0xE0100007L) [static]`

10.54.4.44 WPAError_ScanningIsPending `const uint32_t WPAError_ScanningIsPending = ((0x↵A0220000L) | 0x0036) [static]`

10.54.5 Property Documentation


10.54.5.1 SerialNumber `virtual String^ SerialNumber [get]`

10.55 CMcsUsbPointerContainer Class Reference

10.56 CMEA2100x256FunctionNet Class Reference

[CMEA2100x256FunctionNet](#) is the class to control the MEA2100-256 device needs `#include "Stg200xNet.h"` to resolve documentation reference

Inheritance diagram for CMEA2100x256FunctionNet:

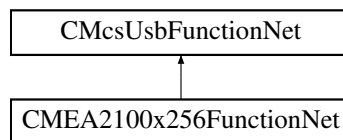

Public Member Functions

- [CMEA2100x256FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pMEA2100x256FunctionPointerContainer)
Initializes a new instance of the [CMEA2100x256FunctionNet](#) class.
- [CMEA2100x256FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMEA2100x256FunctionNet](#) ()
- [ICMEA2100x256FunctionNet](#) ()
- [StimulationLayoutConfigurationEnumNet](#) [GetLayoutConfiguration](#) ()
Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D↔ AC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).
- void [SetLayoutConfiguration](#) ([StimulationLayoutConfigurationEnumNet](#) LayoutConfiguration)
Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D↔ AC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).

Additional Inherited Members

10.56.1 Detailed Description

[CMEA2100x256FunctionNet](#) is the class to control the MEA2100-256 device needs `#include "Stg200xNet.h"` to resolve documentation reference

10.56.2 Constructor & Destructor Documentation

10.56.2.1 CMEA2100x256FunctionNet() [1/2] `CMEA2100x256FunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pMEA2100x256FunctionPointerContainer)`

Initializes a new instance of the [CMEA2100x256FunctionNet](#) class.

10.56.2.2 CMEA2100x256FunctionNet() [2/2] `CMEA2100x256FunctionNet (CMcsUsbNet^ mcsusb)`

10.56.2.3 ~CMEA2100x256FunctionNet() `virtual ~CMEA2100x256FunctionNet () [virtual]`

10.56.2.4 "!CMEA2100x256FunctionNet() `!CMEA2100x256FunctionNet ()`

10.56.3 Member Function Documentation

10.56.3.1 GetLayoutConfiguration() `StimulationLayoutConfigurationEnumNet GetLayoutConfiguration ()`

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D↔ AC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).

Returns

The currently active stimulation layout configuration.

10.56.3.2 SetLayoutConfiguration() `void SetLayoutConfiguration (StimulationLayoutConfigurationEnumNet LayoutConfiguration)`


Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D↔ AC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).

Parameters

<i>LayoutConfiguration</i>	The new stimulation layout configuration.
----------------------------	---

10.57 CMeaAudioFunctionNet Class Reference

Inheritance diagram for CMeaAudioFunctionNet:

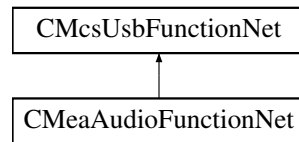

Classes

- struct [s_setaudionet](#)

Public Member Functions

- [CMeaAudioFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meaAudioFunctionPointerContainer)
- [CMeaAudioFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual uint32_t [GetNumberOfAudioChannels](#) ()
Gets the number of available audio channels.
- virtual uint32_t [SetAudioChannels](#) (array< [s_setaudionet](#)[^]> channels)
Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [SetAudioChannels](#) (array< [s_setaudionet](#)[^]> channels, unsigned int virtualDevice)
Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [GetAudioChannels](#) ([System::Runtime::InteropServices::Out]array< [s_setaudionet](#)[^]> channels)
Gets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [GetAudioChannels](#) ([System::Runtime::InteropServices::Out]array< [s_setaudionet](#)[^]> channels, unsigned int virtualDevice)
Gets the electrode to monitor and amplification for the audio channels.

Additional Inherited Members

10.57.1 Constructor & Destructor Documentation

10.57.1.1 CMeaAudioFunctionNet() [1/2] [CMeaAudioFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] meaAudioFunctionPointerContainer)

10.57.1.2 CMeaAudioFunctionNet() [2/2] [CMeaAudioFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.57.2 Member Function Documentation

10.57.2.1 GetAudioChannels() [1/2] virtual uint32_t GetAudioChannels (
 [System::Runtime::InteropServices::Out] array< [s_setaudionet](#)^>^% *channels*)
 [virtual]

Gets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which contains the electrode (channel) and amplification on return.
-----------------	--

Returns

Error Status. 0 on success.

10.57.2.2 GetAudioChannels() [2/2] virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out] array< s_setaudionet^>^% channels, unsigned int virtualDevice) [virtual]

Gets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which contains the electrode (channel) and amplification on return.
-----------------	--

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------

Returns

Error Status. 0 on success.

10.57.2.3 GetNumberOfAudioChannels() virtual uint32_t GetNumberOfAudioChannels () [virtual]

Gets the number of available audio channels.

Returns

The number of audio channels available, 0 when there are none.

10.57.2.4 SetAudioChannels() [1/2] virtual uint32_t SetAudioChannels (array< s_setaudionet^>^ channels) [virtual]

Sets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which defines the electrode (channel) and amplification.
-----------------	---

Returns

Error Status. 0 on success.

10.57.2.5 SetAudioChannels() [2/2] `virtual uint32_t SetAudioChannels (`
 `array< s_setaudionet>^ channels,`
 `unsigned int virtualDevice) [virtual]`

Sets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which defines the electrode (channel) and amplification.
-----------------	---

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------


Returns

Error Status. 0 on success.

10.58 CMeaCleanDeviceNet Class Reference

[CMeaCleanDeviceNet](#) is the class to access the MEA Clean device.

Inheritance diagram for CMeaCleanDeviceNet:

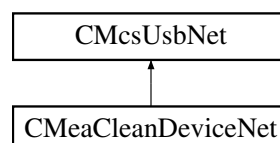

Public Member Functions

- [CMeaCleanDeviceNet](#) ()
Initializes a new instance of the [CMeaCleanDeviceNet](#) class.
- virtual [~CMeaCleanDeviceNet](#) ()
- [!CMeaCleanDeviceNet](#) ()
- void [Start](#) ()
Starts a MEA Clean run.
- void [Stop](#) ()
Stops a MEA Clean run.
- void [SetSlope](#) (uint32_t voltageSlope)
Sets the voltage slope.
- void [SetCycles](#) (uint32_t cycles)
Sets the number of cycles.
- void [SetMinVoltage](#) (int32_t voltageMin)
Sets the lower voltage level.
- void [SetMaxVoltage](#) (int32_t voltageMax)
Sets the upper voltage level.
- bool [IsRunning](#) ()
Gets if the MEA Clean device is running.
- uint32_t [GetSlope](#) ()
Gets the voltage slope.
- uint32_t [GetCycles](#) ()
Gets the number of cycles.
- int32_t [GetMinVoltage](#) ()
Gets the lower voltage level.
- int32_t [GetMaxVoltage](#) ()
Gets the upper voltage level.
- int32_t [GetOutputVoltage](#) ()
Gets the output voltage.
- int32_t [GetCycle](#) ()
Gets the current cycle.

Additional Inherited Members**10.58.1 Detailed Description**

[CMeaCleanDeviceNet](#) is the class to access the MEA Clean device.

10.58.2 Constructor & Destructor Documentation**10.58.2.1 CMeaCleanDeviceNet() [CMeaCleanDeviceNet](#) ()**

Initializes a new instance of the [CMeaCleanDeviceNet](#) class.

10.58.2.2 `~CMeaCleanDeviceNet()` `virtual ~CMeaCleanDeviceNet () [virtual]`

10.58.2.3 `!CMeaCleanDeviceNet()` `!CMeaCleanDeviceNet ()`

10.58.3 Member Function Documentation

10.58.3.1 `GetCycle()` `int32_t GetCycle ()`

Gets the current cycle.

Returns

The cycle number.

10.58.3.2 `GetCycles()` `uint32_t GetCycles ()`

Gets the number of cycles.

Returns

The number of cycles to run for.

10.58.3.3 `GetMaxVoltage()` `int32_t GetMaxVoltage ()`

Gets the upper voltage level

Returns

The upper voltage level in mV.

10.58.3.4 `GetMinVoltage()` `int32_t GetMinVoltage ()`

Gets the lower voltage level.

Returns

The lower voltage level in mV.

10.58.3.5 GetOutputVoltage() `int32_t GetOutputVoltage ()`

Gets the output voltage.

Returns

The output voltage in mV.

10.58.3.6 GetSlope() `uint32_t GetSlope ()`

Gets the voltage slope.

Returns

The voltage slope in mV/s.

10.58.3.7 IsRunning() `bool IsRunning ()`

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

10.58.3.8 SetCycles() `void SetCycles (
uint32_t cycles)`

Sets the number of cycles.

Parameters

<i>cycles</i>	The number of cycles to run for (0 .. 99).
---------------	--

10.58.3.9 SetMaxVoltage() `void SetMaxVoltage (
int32_t voltageMax)`

Sets the upper voltage level.

Parameters

<i>voltageMax</i>	The upper voltage level in mV (-1.6 .. 1.6 V).
-------------------	--

10.58.3.10 SetMinVoltage() `void SetMinVoltage (`
`int32_t voltageMin)`

Sets the lower voltage level.

Parameters

<i>voltageMin</i>	The lower voltage level in mV (-1.6 .. 1.6 V).
-------------------	--

10.58.3.11 SetSlope() `void SetSlope (`
`uint32_t voltageSlope)`

Sets the voltage slope.

Parameters

<i>voltageSlope</i>	The voltage slope in mV/s (range 0 .. 60 V/s).
---------------------	--

10.58.3.12 Start() `void Start ()`

Starts a MEA Clean run.


10.58.3.13 Stop() `void Stop ()`

Stops a MEA Clean run.

10.59 CMeaCoatDeviceNet Class Reference

[CMeaCoatDeviceNet](#) is the class to access the MEA Coat device.

Inheritance diagram for CMeaCoatDeviceNet:

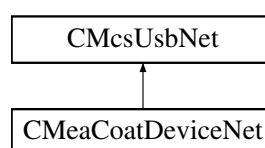

Public Member Functions

- [CMeaCoatDeviceNet](#) ()
Initializes a new instance of the [CMeaCoatDeviceNet](#) class.
- virtual [~CMeaCoatDeviceNet](#) ()
- [!CMeaCoatDeviceNet](#) ()
- void [Start](#) ()
Starts a MEA Coat run.
- void [Stop](#) ()
Stops a MEA Coat run.
- void [SetSlope](#) (int32_t currentSlope)
Sets the current slope.
- void [SetDuration](#) (uint32_t duration)
Sets the duration of a MEA Coat run.
- void [SetMaxCurrent](#) (uint32_t currentMax)
Sets the limit of the current ramp (absolute value).
- void [SetOffsetCurrent](#) (int32_t currentOffset)
Sets the offset of the current.
- bool [IsRunning](#) ()
Gets if the MEA Clean device is running.
- int32_t [GetSlope](#) ()
Gets the current slope.
- uint32_t [GetDuration](#) ()
Gets the duration of a MEA Coat run.
- uint32_t [GetMaxCurrent](#) ()
Gets the limit of the current ramp (absolute value).
- int32_t [GetOffsetCurrent](#) ()
Gets the offset of the current.
- int32_t [GetOutputCurrent](#) ()
Gets the output current.
- int32_t [GetTimeInPlateau](#) ()
Gets the time in the plateau.
- void [SetPauseDuration](#) (uint32_t pauseDuration)
Sets the duration of the pause between MEA Coat pulses.
- uint32_t [GetPauseDuration](#) ()
Gets the duration of the pause between MEA Coat pulses.
- int32_t [GetTimeInPause](#) ()
Gets the time in the pause.
- void [SetCycles](#) (uint32_t cycles)
Sets the number of cycles.
- uint32_t [GetCycles](#) ()
Gets the number of cycles.
- int32_t [GetCurrentCycle](#) ()
Gets the current cycle.

Additional Inherited Members

10.59.1 Detailed Description

[CMeaCoatDeviceNet](#) is the class to access the MEA Coat device.

10.59.2 Constructor & Destructor Documentation

10.59.2.1 CMeaCoatDeviceNet() [CMeaCoatDeviceNet](#) ()

Initializes a new instance of the [CMeaCoatDeviceNet](#) class.

10.59.2.2 ~CMeaCoatDeviceNet() [virtual](#) [~CMeaCoatDeviceNet](#) () [\[virtual\]](#)

10.59.2.3 "!CMeaCoatDeviceNet() [!CMeaCoatDeviceNet](#) ()

10.59.3 Member Function Documentation

10.59.3.1 GetCurrentCycle() [int32_t](#) [GetCurrentCycle](#) ()

Gets the current cycle.

Returns

The cycle number.

10.59.3.2 GetCycles() [uint32_t](#) [GetCycles](#) ()

Gets the number of cycles.

Returns

The number of cycles to run for.

10.59.3.3 GetDuration() [uint32_t](#) [GetDuration](#) ()

Gets the duration of a MEA Coat run.

Returns

The duration in ms.

10.59.3.4 GetMaxCurrent() `uint32_t GetMaxCurrent ()`

Gets the limit of the current ramp (absolute value).

Returns

The limit of the current ramp in pA (absolute value).

10.59.3.5 GetOffsetCurrent() `int32_t GetOffsetCurrent ()`

Gets the offset of the current.

Returns

The offset of the current in pA.

10.59.3.6 GetOutputCurrent() `int32_t GetOutputCurrent ()`

Gets the output current.

Returns

The output current in pA.

10.59.3.7 GetPauseDuration() `uint32_t GetPauseDuration ()`

Gets the duration of the pause between MEA Coat pulses.

Returns

The duration in ms.

10.59.3.8 GetSlope() `int32_t GetSlope ()`

Gets the current slope.

Returns

The current slope in pA/s.

10.59.3.9 GetTimeInPause() `int32_t GetTimeInPause ()`

Gets the time in the pause.

Returns

The time in the pause in ms.

10.59.3.10 GetTimeInPlateau() `int32_t GetTimeInPlateau ()`

Gets the time in the plateau.

Returns

The time in the plateau in ms.

10.59.3.11 IsRunning() `bool IsRunning ()`

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

10.59.3.12 SetCycles() `void SetCycles (
uint32_t cycles)`

Sets the number of cycles.

Parameters

<i>cycles</i>	The number of cycles to run for (0 .. 99).
---------------	--

10.59.3.13 SetDuration() `void SetDuration (
uint32_t duration)`

Sets the duration of a MEA Coat run.

Parameters

<i>duration</i>	The duration in ms (range 0 .. 65 s).
-----------------	---------------------------------------

10.59.3.14 SetMaxCurrent() `void SetMaxCurrent (`
 `uint32_t currentMax)`

Sets the limit of the current ramp (absolute value).

Parameters

<i>currentMax</i>	The limit of the current ramp in pA (absolute value, 0 .. 18 nA).
-------------------	---

10.59.3.15 SetOffsetCurrent() `void SetOffsetCurrent (`
 `int32_t currentOffset)`

Sets the offset of the current.

Parameters

<i>currentOffset</i>	The offset of the current in pA (-10 .. 10 nA).
----------------------	---

10.59.3.16 SetPauseDuration() `void SetPauseDuration (`
 `uint32_t pauseDuration)`

Sets the duration of the pause between MEA Coat pulses.

Parameters

<i>pauseDuration</i>	The duration in ms (range 0 .. 65 s).
----------------------	---------------------------------------

10.59.3.17 SetSlope() `void SetSlope (`
 `int32_t currentSlope)`

Sets the current slope.

Parameters

<i>currentSlope</i>	The current slope in pA/s (range -65 .. 65 nA/s).
---------------------	---

10.59.3.18 Start() `void Start ()`

Starts a MEA Coat run.


10.59.3.19 Stop() void Stop ()

Stops a MEA Coat run.

10.60 CMeaDeviceNet Class Reference

Base class for MEA data acquisition devices.

Inheritance diagram for CMeaDeviceNet:

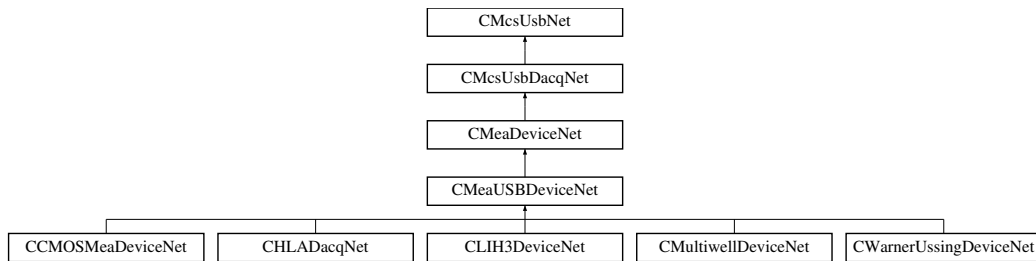

Public Member Functions

- [CMeaDeviceNet](#) (McsBusTypeEnumNet bustype)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [CMeaDeviceNet](#) (McsBusTypeEnumNet bustype, [OnChannelData](#)[^] channelData, [OnError](#)[^] error)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [~CMeaDeviceNet](#) ()
- virtual int32_t [GetGain](#) ()
Gets the amplifier gain of the device.
- int32_t [GetEnumerationSpeed](#) ()
- virtual int32_t [GetAnalogGain](#) ()
Gets the gain of the analog inputs of the device.
- virtual uint32_t [EnableDigitalIn](#) (bool enable, unsigned int virtualDevice)
Enable the digital data word in the datastream.
- virtual uint32_t [EnableDigitalIn](#) (DigitalDatastreamEnableEnumNet enable, unsigned int virtualDevice)
Enable digital data words in the datastream.
- virtual uint32_t [EnableTimestamp](#) (bool enable, unsigned int virtualDevice)
Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.
- virtual uint32_t [EnableChecksum](#) (bool enable, unsigned int virtualDevice)
Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.
- virtual void [SetDigitalOut](#) (unsigned int digout_value, int pulselength)
Generate a pulse on the digital output.
- virtual uint32_t [SetNumberOfChannels](#) (int NumberOfChannels)
Sets the number of analog channels in the datastream.
- virtual uint32_t [SetNumberOfChannels](#) (int NumberOfChannels, unsigned int virtualDevice)
Sets the number of analog channels in the datastream.

- virtual uint32_t [SetNumberOfAnalogChannels](#) (unsigned int NumberOfChannels_HS1, unsigned int NumberOfChannels_HS2, unsigned int NumberOfChannels_DSP, unsigned int NumberOfChannels_IF, unsigned int virtualDevice)
Sets the number of analog channels in the datastream for the MEA2100 device.
- virtual uint32_t [SetTriggerPeriod](#) (int samples, unsigned int virtualDevice)
Sets the maximum number of samples per trigger.
- virtual uint32_t [SetTriggerMaskValue](#) (unsigned int mask, unsigned int value, unsigned int virtualDevice)
Defines a pattern on the digital dataword which will start a trigger when found.

Properties

- [CMeFunctionNet](#)^ [MeFunctionNet](#) [get]
- [CWClassicFunctionNet](#)^ [WClassicFunctionNet](#) [get]
- [CW2100_FunctionNet](#)^ [W2100_FunctionNet](#) [get]
- [CMeaAudioFunctionNet](#)^ [MeaAudioFunctionNet](#) [get]
- [CMeaDigitalDataFunctionNet](#)^ [MeaDigitalDataFunctionNet](#) [get]
- [CMeaFeedbackFunctionNet](#)^ [MeaFeedbackFunctionNet](#) [get]
- virtual int [Gain](#) [get]
The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.
- virtual int [AnalogGain](#) [get]
The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Additional Inherited Members

10.60.1 Detailed Description

Base class for MEA data acquisition devices.

There are two different device types for MEA data acquisition devices. There are the USB-MEA devices and the MC↔_Card. In .NET both classes can be accessed by the constructor of the base class [CMeaDeviceNet](#), which constructs the correct underlying C++ class for the USB-MEA device on the one hand or the MC_Card device on the other hand. Through this interface both device types USB-MEA devices and MC_Card devices can be accessed

10.60.2 Constructor & Destructor Documentation

10.60.2.1 [CMeaDeviceNet\(\)](#) [1/2] [CMeaDeviceNet](#) (McsBusTypeEnumNet *bustype*)

Initializes a new instance of [CMeaDeviceNet](#) class.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

10.60.2.2 CMeaDeviceNet() [2/2] [CMeaDeviceNet](#) (
 McsBusTypeEnumNet *bustype*,
 [OnChannelData](#)[^] *channelData*,
 [OnError](#)[^] *error*)

Initializes a new instance of [CMeaDeviceNet](#) class.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

Parameters

<i>channelData</i>	Callback to call when new data is available.
--------------------	--

Parameters

<i>error</i>	Callback to call when an error occurred.
--------------	--

10.60.2.3 ~CMeaDeviceNet() [~CMeaDeviceNet](#) ()

10.60.3 Member Function Documentation

10.60.3.1 EnableChecksum() virtual uint32_t EnableChecksum (
 bool *enable*,
 unsigned int *virtualDevice*) [virtual]

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

10.60.3.2 EnableDigitalIn() [1/2] virtual uint32_t EnableDigitalIn (
 bool *enable*,
 unsigned int *virtualDevice*) [virtual]

Enable the digital data word in the datastream.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

10.60.3.3 EnableDigitalIn() [2/2] virtual uint32_t EnableDigitalIn (
 DigitalDatastreamEnableEnumNet *enable*,
 unsigned int *virtualDevice*) [virtual]

Enable digital data words in the datastream.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

10.60.3.4 EnableTimestamp() virtual uint32_t EnableTimestamp (
 bool *enable*,
 unsigned int *virtualDevice*) [virtual]

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

10.60.3.5 GetAnalogGain() `virtual int32_t GetAnalogGain () [virtual]`

Gets the gain of the analog inputs of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

10.60.3.6 GetEnumerationSpeed() `int32_t GetEnumerationSpeed ()`**10.60.3.7 GetGain()** `virtual int32_t GetGain () [virtual]`

Gets the amplifier gain of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

10.60.3.8 SetDigitalOut() `virtual void SetDigitalOut (
 unsigned int digout_value,
 int pulselength) [virtual]`

Generate a pulse on the digital output.

Parameters

<i>digout_value</i>	Bitmask to set on the digital out.
---------------------	------------------------------------

Parameters

<i>pulselength</i>	Pulselength in ms.
--------------------	--------------------

10.60.3.9 SetNumberOfAnalogChannels() virtual uint32_t SetNumberOfAnalogChannels (
 unsigned int *NumberOfChannels_HS1*,
 unsigned int *NumberOfChannels_HS2*,
 unsigned int *NumberOfChannels_DSP*,
 unsigned int *NumberOfChannels_IF*,
 unsigned int *virtualDevice*) [virtual]

Sets the number of analog channels in the datastream for the MEA2100 device.

Parameters

<i>NumberOfChannels_HS1</i>	Number of analog channels from the Headstage 1.
-----------------------------	---

Parameters

<i>NumberOfChannels_HS2</i>	Number of analog channels from the Headstage 2.
-----------------------------	---

Parameters

<i>NumberOfChannels_DSP</i>	Number of data words from the DSP.
-----------------------------	------------------------------------

Parameters

<i>NumberOfChannels_IF</i>	Number of analog channels from the Interfaceboard.
----------------------------	--

Parameters

<i>virtualDevice</i>	virtualDevice to use.
----------------------	-----------------------

Returns

Error Status. 0 on success.

10.60.3.10 SetNumberOfChannels() [1/2] virtual uint32_t SetNumberOfChannels (
int *NumberOfChannels*) [virtual]

Sets the number of analog channels in the datastream.

Parameters

<i>NumberOfChannels</i>	Number of analog channels.
-------------------------	----------------------------

Returns

Error Status. 0 on success.

10.60.3.11 SetNumberOfChannels() [2/2] virtual uint32_t SetNumberOfChannels (
int *NumberOfChannels*,
unsigned int *virtualDevice*) [virtual]

Sets the number of analog channels in the datastream.

Parameters

<i>NumberOfChannels</i>	Number of analog channels.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

10.60.3.12 SetTriggerMaskValue() virtual uint32_t SetTriggerMaskValue (
unsigned int *mask*,
unsigned int *value*,
unsigned int *virtualDevice*) [virtual]

Defines a pattern on the digital dataword which will start a trigger when found.

Parameters

<i>mask</i>	Bits in the digital dataword which are monitored for a match with value.
-------------	--

Parameters

<i>value</i>	Pattern which must match for the trigger to start.
--------------	--

Returns

Error Status. 0 on success.

10.60.3.13 SetTriggerPeriod() `virtual uint32_t SetTriggerPeriod (`
`int samples,`
`unsigned int virtualDevice) [virtual]`

Sets the maximum number of samples per trigger.

Parameters

<i>samples</i>	Number of samples to acquire after the trigger condition is met.
----------------	--

Returns

Error Status. 0 on success.

10.60.4 Property Documentation

10.60.4.1 AnalogGain `virtual int AnalogGain [get]`

The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

10.60.4.2 Gain `virtual int Gain [get]`

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

10.60.4.3 MeaAudioFunctionNet `CMeaAudioFunctionNet^ MeaAudioFunctionNet [get]`

10.60.4.4 MeaDigitalDataFunctionNet `CMeaDigitalDataFunctionNet^ MeaDigitalDataFunctionNet [get]`

10.60.4.5 MeaFeedbackFunctionNet [CMeaFeedbackFunctionNet](#)[^] MeaFeedbackFunctionNet [get]


10.60.4.6 MeFunctionNet [CMeFunctionNet](#)[^] MeFunctionNet [get]

10.60.4.7 W2100_FunctionNet [CW2100_FunctionNet](#)[^] W2100_FunctionNet [get]

10.60.4.8 WClassicFunctionNet [CWClassicFunctionNet](#)[^] WClassicFunctionNet [get]

10.61 CMeaDigitalDataFunctionNet Class Reference

Inheritance diagram for CMeaDigitalDataFunctionNet:

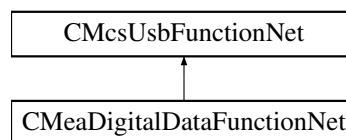

Public Member Functions

- [CMeaDigitalDataFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meaDigitalFunctionPointerContainer)
- [CMeaDigitalDataFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetDigitalData](#) (unsigned int digital_value, unsigned int digital_value_mask)
Generate a value on the digital output.
- void [SetDigitalData](#) (unsigned int bit_number, bool value)
Generate a value on the digital output.
- unsigned int [GetDigitalData](#) ()
Get the value of the digital output.

Additional Inherited Members

10.61.1 Constructor & Destructor Documentation

10.61.1.1 CMeaDigitalDataFunctionNet() [1/2] [CMeaDigitalDataFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] meaDigitalFunctionPointerContainer)

10.61.1.2 CMeaDigitalDataFunctionNet() [2/2] `CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb)`

10.61.2 Member Function Documentation

10.61.2.1 GetDigitalData() `unsigned int GetDigitalData ()`

Get the value of the digital output.

Returns

Value on the digital data register.

10.61.2.2 SetDigitalData() [1/2] `void SetDigitalData (unsigned int bit_number, bool value)`

Generate a value on the digital output.

Parameters

<i>bit_number</i>	Bit number to change.
-------------------	-----------------------

Parameters

<i>value</i>	Bit value.
--------------	------------

10.61.2.3 SetDigitalData() [2/2] `void SetDigitalData (unsigned int digital_value, unsigned int digital_value_mask)`

Generate a value on the digital output.

Parameters


<i>digital_value</i>	Value to set.
----------------------	---------------

Parameters

<i>digital_value_mask</i>	Mask for change.
---------------------------	------------------

10.62 CMeaFeedbackFunctionNet Class Reference

Inheritance diagram for CMeaFeedbackFunctionNet:

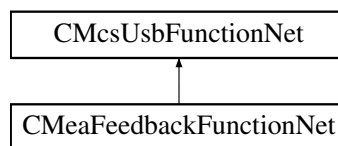

Public Member Functions

- [CMeaFeedbackFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meaFeedback↔FunctionNet)
- [CMeaFeedbackFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [FeedbackSetFeedback](#) (unsigned char on, unsigned short digoutmask, unsigned short diginmask)
- unsigned int [FeedbackGetSampleTimerCount](#) ([System::Runtime::InteropServices::Out]unsigned int% CurrentCount, [System::Runtime::InteropServices::Out]unsigned int% LastKnownCount, [System::Runtime::InteropServices::Out]bool% On)
- void [FeedbackSetDigitalMapping](#) (unsigned short channel, unsigned short outmapping, unsigned short inmapping)
- void [FeedbackSetFilterParameter](#) (unsigned char filter, array< short >[^] parameters)
- void [FeedbackSetFilterParameter32](#) (unsigned char filter, array< int >[^] parameters)
- void [FeedbackSetIIRFilterParameter](#) (unsigned char filter, int length, array< double >[^] parameters)
- void [FeedbackSetMkFilter](#) (unsigned char filter, String[^] filtertype, double cheb_ripple, String[^] passtype, int order, double alpha1, double alpha2)
- void [FeedbackSetChannelFilter](#) (short channel, char filter)
- void [FeedbackSetGlobalChannelFilter](#) (char filter, unsigned short firstchannel, unsigned short lastchannel)
- void [FeedbackSetFilterOff](#) ()
- void [FeedbackSetNumberOfSpikeDetectors](#) (unsigned short number)
- void [FeedbackSetSpikeDetectorThreshold](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, int threshold1, int threshold2, short slope)
- void [FeedbackSetNumberOfRateCounter](#) (unsigned short number)
- void [FeedbackSetRateCounter](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel)
- void [FeedbackSetNumberOfRateDetectors](#) (unsigned short number)
- void [FeedbackSetRateDetector](#) (unsigned short position, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, unsigned short pulses, unsigned int duration1, unsigned int duration2)
- void [FeedbackSetNumberOfLogics](#) (unsigned short number)
- void [FeedbackSetLogic](#) (unsigned short position, array< unsigned short >[^] sourcechannel, unsigned short resultchannel, unsigned int lookup)
- void [FeedbackSetNumberOfTriggers](#) (unsigned short number)
- void [FeedbackSetTrigger](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit)
- void [FeedbackSetAnalogSource](#) (AnalogSourceEnumNet AnalogSource, unsigned int Channels, unsigned int Offset)

Additional Inherited Members

10.62.1 Constructor & Destructor Documentation

10.62.1.1 CMeaFeedbackFunctionNet() [1/2] `CMeaFeedbackFunctionNet` (
`CMcsUsbNet`[^] *mcsusb*,
`CMcsUsbFunctionPointerContainer`[^] *meaFeedbackFunctionNet*)

10.62.1.2 CMeaFeedbackFunctionNet() [2/2] `CMeaFeedbackFunctionNet` (
`CMcsUsbNet`[^] *mcsusb*)

10.62.2 Member Function Documentation

10.62.2.1 FeedbackGetSampleTimerCount() `unsigned int FeedbackGetSampleTimerCount` (
`[System::Runtime::InteropServices::Out]` `unsigned int` % *CurrentCount*,
`[System::Runtime::InteropServices::Out]` `unsigned int` % *LastKnownCount*,
`[System::Runtime::InteropServices::Out]` `bool` % *On*)

10.62.2.2 FeedbackSetAnalogSource() `void FeedbackSetAnalogSource` (
`AnalogSourceEnumNet` *AnalogSource*,
`unsigned int` *Channels*,
`unsigned int` *Offset*)

10.62.2.3 FeedbackSetChannelFilter() `void FeedbackSetChannelFilter` (
`short` *channel*,
`char` *filter*)

10.62.2.4 FeedbackSetDigitalMapping() `void FeedbackSetDigitalMapping` (
`unsigned short` *channel*,
`unsigned short` *outmapping*,
`unsigned short` *inmapping*)

10.62.2.5 FeedbackSetFeedback() void FeedbackSetFeedback (
 unsigned char *on*,
 unsigned short *digoutmask*,
 unsigned short *diginmask*)

10.62.2.6 FeedbackSetFilterOff() void FeedbackSetFilterOff ()

10.62.2.7 FeedbackSetFilterParameter() void FeedbackSetFilterParameter (
 unsigned char *filter*,
 array< short >^ *parameters*)

10.62.2.8 FeedbackSetFilterParameter32() void FeedbackSetFilterParameter32 (
 unsigned char *filter*,
 array< int >^ *parameters*)

10.62.2.9 FeedbackSetGlobalChannelFilter() void FeedbackSetGlobalChannelFilter (
 char *filter*,
 unsigned short *firstchannel*,
 unsigned short *lastchannel*)

10.62.2.10 FeedbackSetIIRFilterParameter() void FeedbackSetIIRFilterParameter (
 unsigned char *filter*,
 int *length*,
 array< double >^ *parameters*)

10.62.2.11 FeedbackSetLogic() void FeedbackSetLogic (
 unsigned short *position*,
 array< unsigned short >^ *sourcechannel*,
 unsigned short *resultchannel*,
 unsigned int *lookup*)

10.62.2.12 FeedbackSetMkFilter() void FeedbackSetMkFilter (
 unsigned char *filter*,
 String^ *filtertype*,
 double *cheb_ribble*,
 String^ *passtype*,
 int *order*,
 double *alpha1*,
 double *alpha2*)

10.62.2.13 FeedbackSetNumberOfLogics() void FeedbackSetNumberOfLogics (
 unsigned short *number*)

10.62.2.14 FeedbackSetNumberOfRateCounter() void FeedbackSetNumberOfRateCounter (
 unsigned short *number*)

10.62.2.15 FeedbackSetNumberOfRateDetectors() void FeedbackSetNumberOfRateDetectors (
 unsigned short *number*)

10.62.2.16 FeedbackSetNumberOfSpikeDetectors() void FeedbackSetNumberOfSpikeDetectors (
 unsigned short *number*)

10.62.2.17 FeedbackSetNumberOfTriggers() void FeedbackSetNumberOfTriggers (
 unsigned short *number*)

10.62.2.18 FeedbackSetRateCounter() void FeedbackSetRateCounter (
 unsigned short *position*,
 unsigned short *sourcechannel*,
 unsigned short *resultchannel*)


10.62.2.19 FeedbackSetRateDetector() void FeedbackSetRateDetector (
 unsigned short *position*,
 unsigned short *resultchannel*,
 unsigned short *trigger*,
 unsigned short *totzeit*,
 unsigned short *pulses*,
 unsigned int *duration1*,
 unsigned int *duration2*)

10.62.2.20 FeedbackSetSpikeDetectorThreshold() `void FeedbackSetSpikeDetectorThreshold (`
 `unsigned short position,`
 `unsigned short sourcechannel,`
 `unsigned short resultchannel,`
 `unsigned short trigger,`
 `unsigned short totzeit,`
 `int threshold1,`
 `int threshold2,`
 `short slope)`

10.62.2.21 FeedbackSetTrigger() `void FeedbackSetTrigger (`
 `unsigned short position,`
 `unsigned short sourcechannel,`
 `unsigned short resultchannel,`
 `unsigned short trigger,`
 `unsigned short totzeit)`

10.63 CMealImpedanceDeviceNet Class Reference

Inheritance diagram for CMealImpedanceDeviceNet:

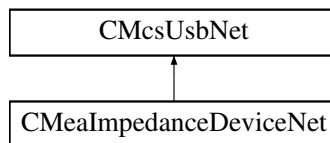

Public Member Functions

- [CMealImpedanceDeviceNet](#) ()
- [~CMealImpedanceDeviceNet](#) ()
- void [StartMeasurement](#) (unsigned short channel)
- unsigned short [GetReady](#) ()
- unsigned short [GetArraySize](#) ()
- array< unsigned short > ^ [GetResult](#) ()
- unsigned short [GetAdapterCode](#) ()
- unsigned int [GetImpedanceTestFrequency](#) ()
- void [SetImpedanceTestFrequency](#) (unsigned int TestFrequency_Hertz)

Additional Inherited Members

10.63.1 Constructor & Destructor Documentation

10.63.1.1 CMealImpedanceDeviceNet() [CMealImpedanceDeviceNet](#) ()

10.63.1.2 `~CMeaImpedanceDeviceNet()` `~CMeaImpedanceDeviceNet ()`

10.63.2 Member Function Documentation

10.63.2.1 `GetAdapterCode()` `unsigned short GetAdapterCode ()`

10.63.2.2 `GetArraySize()` `unsigned short GetArraySize ()`

10.63.2.3 `GetImpedanceTestFrequency()` `unsigned int GetImpedanceTestFrequency ()`

10.63.2.4 `GetReady()` `unsigned short GetReady ()`

10.63.2.5 `GetResult()` `array<unsigned short> ^ GetResult ()`


10.63.2.6 `SetImpedanceTestFrequency()` `void SetImpedanceTestFrequency (`
`unsigned int TestFrequency_Hertz)`

10.63.2.7 `StartMeasurement()` `void StartMeasurement (`
`unsigned short channel)`

10.64 CMeasureTableDeviceNet Class Reference

[CMeasureTableDeviceNet](#) is the to control the MCS HLA device

Inheritance diagram for CMeasureTableDeviceNet:

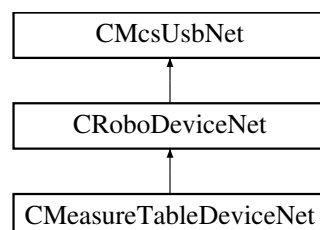

Public Member Functions

- [CMeasureTableDeviceNet](#) (void)

Properties

- [CMcsBus_SensorNet^ Sensor](#) [get]

Additional Inherited Members

10.64.1 Detailed Description

[CMeasureTableDeviceNet](#) is the to control the MCS HLA device

10.64.2 Constructor & Destructor Documentation

10.64.2.1 CMeasureTableDeviceNet() [CMeasureTableDeviceNet](#) (
void)


10.64.3 Property Documentation

10.64.3.1 Sensor [CMcsBus_SensorNet^ Sensor](#) [get]

10.65 CMeaSwitchDeviceNet Class Reference

The class to control the USB-MEA-Switch.

Inheritance diagram for CMeaSwitchDeviceNet:

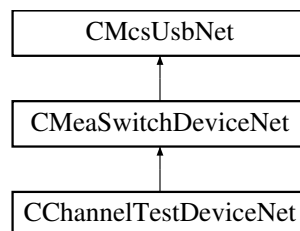

Public Member Functions

- [CMeaSwitchDeviceNet](#) ()
Constructor.
- [~CMeaSwitchDeviceNet](#) ()
Destructor.
- unsigned short [GetNumber](#) ()
Gets the number of boards in the device.
- array< unsigned char > ^ [GetPattern](#) ()
Gets the pattern of the switches that are currently set in the device as char array.
- array< bool > ^ [GetPatternBool](#) ()
Gets the pattern of the switches that are currently set in the device as bools.
- void [SetPattern](#) (array< unsigned char >^ pattern)
Sets the pattern of switches from a char array.
- void [SetPatternBool](#) (array< bool >^ pattern)
Sets the pattern of switches from a.

Additional Inherited Members**10.65.1 Detailed Description**

The class to control the USB-MEA-Switch.

This class controls the settings of the USB-MEA-Switch. The box has two inputs for signals from a MEA amplifier. Each of the 64 outputs can be connected to one of the MEAs at the same channel.

10.65.2 Constructor & Destructor Documentation**10.65.2.1 CMeaSwitchDeviceNet()** [CMeaSwitchDeviceNet](#) ()

Constructor.

10.65.2.2 ~CMeaSwitchDeviceNet() [~CMeaSwitchDeviceNet](#) ()

Destructor.

10.65.3 Member Function Documentation**10.65.3.1 GetNumber()** unsigned short [GetNumber](#) ()

Gets the number of boards in the device.

The MEA-Switch are delivered with 64 or 128 channels

10.65.3.2 GetPattern() `array<unsigned char> ^ GetPattern ()`

Gets the pattern of the switches that are currently set in the device as char array.

10.65.3.3 GetPatternBool() `array<bool> ^ GetPatternBool ()`

Gets the pattern of the switches that are currently set in the device as bools.

10.65.3.4 SetPattern() `void SetPattern (array< unsigned char > ^ pattern)`

Sets the pattern of switches from a char array.


10.65.3.5 SetPatternBool() `void SetPatternBool (array< bool > ^ pattern)`

Sets the pattern of switches from a.

10.66 CMeaUSBDeviceNet Class Reference

Class for data acquisition via ME and MEA USB amplifiers

Inheritance diagram for CMeaUSBDeviceNet:

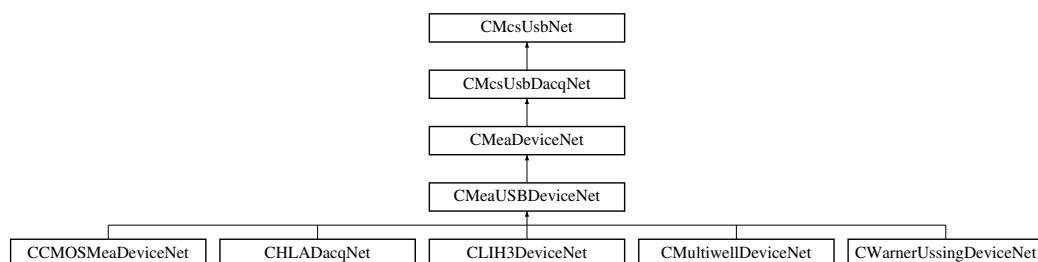

Public Member Functions

- [CMeaUSBDeviceNet](#) ([OnChannelData](#) ^ channelData, [OnError](#) ^ error)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [CMeaUSBDeviceNet](#) ()
Initializes a new instance of [CMeaDeviceNet](#) class.
- [~CMeaUSBDeviceNet](#) ()

Additional Inherited Members

10.66.1 Detailed Description

Class for data acquisition via ME and MEA USB amplifiers

10.66.2 Constructor & Destructor Documentation

10.66.2.1 CMeaUSBDeviceNet() [1/2] CMeaUSBDeviceNet (OnChannelData^ *channelData*, OnError^ *error*)

Initializes a new instance of [CMeaDeviceNet](#) class.

Parameters

<i>channelData</i>	Handler to call when new data is available.
--------------------	---

Parameters

<i>error</i>	Handler to call when an error occurs.
--------------	---------------------------------------


10.66.2.2 CMeaUSBDeviceNet() [2/2] CMeaUSBDeviceNet ()

Initializes a new instance of [CMeaDeviceNet](#) class.

10.66.2.3 ~CMeaUSBDeviceNet() ~CMeaUSBDeviceNet ()

10.67 CMeFunctionNet Class Reference

Inheritance diagram for CMeFunctionNet:

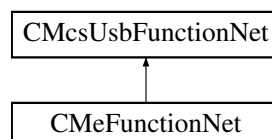

Public Member Functions

- [CMeFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meFunctionPointer↔ Container)
Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.
- [CMeFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMeFunctionNet](#) (void)
- [!CMeFunctionNet](#) (void)
- void [SetTransformer](#) (unsigned int index, bool onoff)

Additional Inherited Members

10.67.1 Detailed Description

10.67.2 Constructor & Destructor Documentation

10.67.2.1 CMeFunctionNet() [1/2] [CMeFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] meFunctionPointerContainer)

Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.

10.67.2.2 CMeFunctionNet() [2/2] [CMeFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.67.2.3 ~CMeFunctionNet() virtual [~CMeFunctionNet](#) (
void) [virtual]

10.67.2.4 !CMeFunctionNet() [!CMeFunctionNet](#) (
void)


10.67.3 Member Function Documentation

10.67.3.1 SetTransformer() void [SetTransformer](#) (
unsigned int index,
bool onoff)

10.68 CMultiBatteryChargerDeviceNet Class Reference

[CMultiBatteryChargerDeviceNet](#) is the class to access the MBC-08 device.

Inheritance diagram for CMultiBatteryChargerDeviceNet:

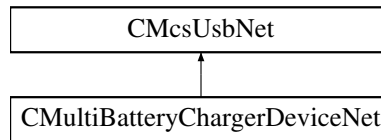

Public Member Functions

- [CMultiBatteryChargerDeviceNet](#) ()
Initializes a new instance of the [CMultiBatteryChargerDeviceNet](#) class.
- virtual [~CMultiBatteryChargerDeviceNet](#) ()
- [!CMultiBatteryChargerDeviceNet](#) ()
- [uint32_t GetChargeCurrent](#) (uint32_t NrChannel)
gets the charge current; unit: mA
- [uint32_t GetDischargeCurrent](#) (uint32_t NrChannel)
gets the discharge current; unit: mA
- void [SetDischargeCurrentSetPoint](#) (uint32_t NrChannel, uint32_t DischargeCurrent_mA)
sets the setpoint for the discharge current; unit: mA
- [uint32_t GetDischargeCurrentSetPoint](#) (uint32_t NrChannel)
gets the setpoint for the discharge current; unit: mA
- void [SetFinalDischargeVoltage](#) (uint32_t NrChannel, uint32_t FinalDischargeVoltage_mV)
sets the final discharge voltage; unit: mV
- [uint32_t GetFinalDischargeVoltage](#) (uint32_t NrChannel)
gets the final discharge voltage; unit: mV
- [uint32_t GetDischargeCapacity](#) (uint32_t NrChannel)
gets the discharge capacity; unit: ?Ah
- [uint32_t GetChargeCapacity](#) (uint32_t NrChannel)
gets the charge capacity; unit: ?Ah
- [uint32_t GetBatteryVoltage](#) (uint32_t NrChannel)
gets the battery voltage; unit: mV
- [uint32_t GetChannels](#) ()
gets number of channels
- void [SetRatedCapacityVolatile](#) (uint32_t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)
sets the rated capacity (i.e. charge current) without storing it persistently
- void [SetChargingMode](#) (uint32_t NrChannel, MbcChargingModeEnumNet NewOperatingMode)
sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge
- MbcChargingModeEnumNet [GetChargingMode](#) (uint32_t NrChannel)
gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge
- MbcChannelStateEnumNet [GetChannelState](#) (uint32_t NrChannel)
gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge
- void [CapacityTest](#) (uint32_t NrChannel)
start capacity test on channel
- void [ChannelReset](#) (uint32_t NrChannel)
cancel charging and capacity test functions; check if battery is connected

- void [SetChargingPCoefficient](#) (uint32_t pCoefficient)
sets the p-coefficient for charging in mA/V / nominal charging current
- uint32_t [GetChargingPCoefficient](#) ()
gets the p-coefficient for charging in mA/V / nominal charging current
- void [SetRatedCapacity](#) (uint32_t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)
sets the rated capacity
- MbcRatedCapacityEnumNet [GetRatedCapacity](#) (uint32_t NrChannel)
gets the rated capacity

Additional Inherited Members

10.68.1 Detailed Description

[CMultiBatteryChargerDeviceNet](#) is the class to access the MBC-08 device.

10.68.2 Constructor & Destructor Documentation

10.68.2.1 CMultiBatteryChargerDeviceNet() [CMultiBatteryChargerDeviceNet](#) ()

Initializes a new instance of the [CMultiBatteryChargerDeviceNet](#) class.

10.68.2.2 ~CMultiBatteryChargerDeviceNet() [virtual ~CMultiBatteryChargerDeviceNet](#) () [virtual]

10.68.2.3 "!CMultiBatteryChargerDeviceNet() [!CMultiBatteryChargerDeviceNet](#) ()

10.68.3 Member Function Documentation

10.68.3.1 CapacityTest() [void CapacityTest](#) (uint32_t NrChannel)

start capacity test on channel

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

10.68.3.2 ChannelReset() `void ChannelReset (`
 `uint32_t NrChannel)`

cancel charging and capacity test functions; check if battery is connected

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

10.68.3.3 GetBatteryVoltage() `uint32_t GetBatteryVoltage (`
 `uint32_t NrChannel)`

gets the battery voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the battery voltage in mV

10.68.3.4 GetChannels() `uint32_t GetChannels ()`

gets number of channels

Returns

number of channels

10.68.3.5 GetChannelState() `MbcChannelStateEnumNet GetChannelState (`
 `uint32_t NrChannel)`

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the current state

10.68.3.6 GetChargeCapacity() `uint32_t GetChargeCapacity (`
`uint32_t NrChannel)`

gets the charge capacity; unit: ?Ah

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity in uAh

10.68.3.7 GetChargeCurrent() `uint32_t GetChargeCurrent (`
`uint32_t NrChannel)`

gets the charge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the measured charge current in mA

10.68.3.8 GetChargingMode() `MbcChargingModeEnumNet GetChargingMode (`
`uint32_t NrChannel)`

gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the charging mode

10.68.3.9 GetChargingPCoefficient() `uint32_t GetChargingPCoefficient ()`

gets the p-coefficient for charging in mA/V / nominal charging current

Returns

the p-coefficient

10.68.3.10 GetDischargeCapacity() `uint32_t GetDischargeCapacity (
uint32_t NrChannel)`

gets the discharge capacity; unit: ?Ah

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity in uAh

10.68.3.11 GetDischargeCurrent() `uint32_t GetDischargeCurrent (
uint32_t NrChannel)`

gets the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the measured discharge current in mA

10.68.3.12 GetDischargeCurrentSetPoint() `uint32_t GetDischargeCurrentSetPoint (
uint32_t NrChannel)`

gets the setpoint for the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the discharge current in mA

10.68.3.13 GetFinalDischargeVoltage() `uint32_t GetFinalDischargeVoltage (uint32_t NrChannel)`

gets the final discharge voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the battery voltage in mV at the end of discharge

10.68.3.14 GetRatedCapacity() `MbcRatedCapacityEnumNet GetRatedCapacity (uint32_t NrChannel)`

gets the rated capacity

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity

10.68.3.15 SetChargingMode() `void SetChargingMode (uint32_t NrChannel, MbcChargingModeEnumNet NewOperatingMode)`

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
<i>NewOperatingMode</i>	the charging mode

10.68.3.16 SetChargingPCoefficient() `void SetChargingPCoefficient (`
 `uint32_t pCoefficient)`

sets the p-coefficient for charging in mA/V / nominal charging current

Parameters

<i>pCoefficient</i>	the p-coefficient
---------------------	-------------------

10.68.3.17 SetDischargeCurrentSetPoint() `void SetDischargeCurrentSetPoint (`
 `uint32_t NrChannel,`
 `uint32_t DischargeCurrent_mA)`

sets the setpoint for the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
<i>DischargeCurrent_mA</i>	the discharge current in mA

10.68.3.18 SetFinalDischargeVoltage() `void SetFinalDischargeVoltage (`
 `uint32_t NrChannel,`
 `uint32_t FinalDischargeVoltage_mV)`

sets the final discharge voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
<i>FinalDischargeVoltage_mV</i>	the battery voltage in mV at the end of discharge

10.68.3.19 SetRatedCapacity() `void SetRatedCapacity (`
 `uint32_t NrChannel,`
 `MbcRatedCapacityEnumNet NewRatedCapacity)`

sets the rated capacity

Parameters

<i>NrChannel</i>	the channel number
<i>NewRatedCapacity</i>	the capacity

10.68.3.20 SetRatedCapacityVolatile() void SetRatedCapacityVolatile (
 uint32_t NrChannel,
 MbcRatedCapacityEnumNet NewRatedCapacity)

sets the rated capacity (i.e. charge current) without storing it persistently


Parameters

<i>NrChannel</i>	the channel number
<i>NewRatedCapacity</i>	the capacity

10.69 CMultiwellCallbackFunctionNet Class Reference

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

Inheritance diagram for CMultiwellCallbackFunctionNet:

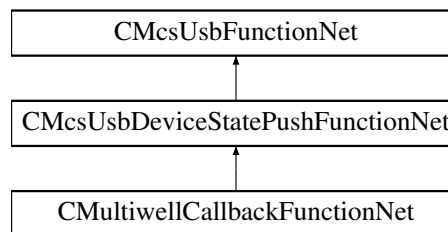

Public Member Functions

- delegate void [OnGetPlateClampStateByHeadstage](#) (uint32_t Headstage, PlateClampEnumNet plateState)
- delegate void [OnGetPlateTypeByHeadstage](#) (uint32_t Headstage, MultiwellPlateTypeEnumNet plateType)
- delegate void [OnIsPlateTypeValidByHeadstage](#) (uint32_t Headstage, bool isValid)
- [CMultiwellCallbackFunctionNet](#) (CMcsUsbNet[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pMultiwell↔ [CallbackFunctionPointerContainer](#))
Initializes a new instance of the [CMultiwellCallbackFunctionNet](#) class.
- [CMultiwellCallbackFunctionNet](#) (CMcsUsbNet[^] mcsusb)
- virtual [~CMultiwellCallbackFunctionNet](#) ()
- [!CMultiwellCallbackFunctionNet](#) ()
- PlateClampEnumNet [GetPlateClampStateByHeadstage](#) (uint32_t Headstage)
Gets the state of the plate
- MultiwellPlateTypeEnumNet [GetPlateTypeByHeadstage](#) (uint32_t Headstage)
Gets the plate type.
- bool [IsPlateTypeValidByHeadstage](#) (uint32_t Headstage)
Checks whether the plate type is valid, meaning all pins have contact.

Events

- [OnGetPlateClampStateByHeadstage](#)[^] [GetPlateClampStateByHeadstageEvent](#) [add, remove, raise]
Event fires when the plate state for the headstage number has changed
- [OnGetPlateTypeByHeadstage](#)[^] [GetPlateTypeByHeadstageEvent](#) [add, remove, raise]
Event fires when the plate type for the headstage to query has changed
- [OnIsPlateTypeValidByHeadstage](#)[^] [IsPlateTypeValidByHeadstageEvent](#) [add, remove, raise]
Event fires when "true" when all pins have contact, otherwise "false" for the headstage to query has changed

Additional Inherited Members

10.69.1 Detailed Description

[CMultiwellCallbackFunctionNet](#) is the class to access the Multiwell-Mini-Stimulator

10.69.2 Constructor & Destructor Documentation

10.69.2.1 CMultiwellCallbackFunctionNet() [1/2] [CMultiwellCallbackFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pMultiwellCallbackFunctionPointerContainer*)

Initializes a new instance of the [CMultiwellCallbackFunctionNet](#) class.

10.69.2.2 CMultiwellCallbackFunctionNet() [2/2] [CMultiwellCallbackFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

10.69.2.3 ~CMultiwellCallbackFunctionNet() `virtual ~CMultiwellCallbackFunctionNet ()` [virtual]

10.69.2.4 !CMultiwellCallbackFunctionNet() `!CMultiwellCallbackFunctionNet ()`

10.69.3 Member Function Documentation

10.69.3.1 GetPlateClampStateByHeadstage() `PlateClampEnumNet GetPlateClampStateByHeadstage (
 uint32_t Headstage)`

Gets the state of the plate

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The plate state

10.69.3.2 GetPlateTypeByHeadstage() `MultiwellPlateTypeEnumNet GetPlateTypeByHeadstage (uint32_t Headstage)`

Gets the plate type.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

the plate type

10.69.3.3 IsPlateTypeValidByHeadstage() `bool IsPlateTypeValidByHeadstage (uint32_t Headstage)`

Checks whether the plate type is valid, meaning all pins have contact.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

"true" when all pins have contact, otherwise "false".

10.69.3.4 OnGetPlateClampStateByHeadstage() `delegate void OnGetPlateClampStateByHeadstage (uint32_t Headstage, PlateClampEnumNet plateState)`

10.69.3.5 OnGetPlateTypeByHeadstage() `delegate void OnGetPlateTypeByHeadstage (uint32_t Headstage, MultiwellPlateTypeEnumNet plateType)`

10.69.3.6 OnIsPlateTypeValidByHeadstage() `delegate void OnIsPlateTypeValidByHeadstage (uint32_t Headstage, bool isValid)`

10.69.4 Event Documentation

10.69.4.1 GetPlateClampStateByHeadstageEvent [OnGetPlateClampStateByHeadstage](#)[^] [GetPlateClampStateByHeadstageEvent](#) [add], [remove], [raise]

Event fires when the plate state for the headstage number has changed

10.69.4.2 GetPlateTypeByHeadstageEvent [OnGetPlateTypeByHeadstage](#)[^] [GetPlateTypeByHeadstageEvent](#) [add], [remove], [raise]

Event fires when the plate type for the headstage to query has changed


10.69.4.3 IsPlateTypeValidByHeadstageEvent [OnIsPlateTypeValidByHeadstage](#)[^] [IsPlateTypeValidByHeadstageEvent](#) [add], [remove], [raise]

Event fires when "true" when all pins have contact, otherwise "false" for the headstage to query has changed

10.70 CMultiwellDeviceNet Class Reference

[CMultiwellDeviceNet](#) is the class to access the HEKA LIH3 device.

Inheritance diagram for CMultiwellDeviceNet:

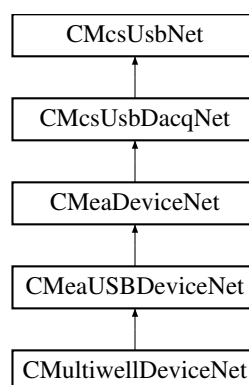

Public Member Functions

- [CMultiwellDeviceNet](#) ()
Initializes a new instance of the [CMultiwellDeviceNet](#) class.
- virtual [~CMultiwellDeviceNet](#) ()
- [ICMultiwellDeviceNet](#) ()
- PlateClampEnumNet [GetPlateClampState](#) ()
Gets the state of the Multiwell plate clamp.
- PlateClampEnumNet [GetPlateClampStateByHeadstage](#) (uint32_t Headstage)
Gets the state of the plate
- void [OpenPlateClamp](#) ()
Opens the plate clamp.
- void [ClosePlateClamp](#) ()
Closes the plate clamp.
- void [StopPlateClamp](#) ()
Stops the plate clamp movement.
- uint32_t [GetPlateClampLockState](#) ()
Gets the state of the plate clamp lock.
- void [LockPlateClamp](#) ()
Locks the plate clamp.
- void [UnlockPlateClamp](#) ()
Unlocks the plate clamp.
- MultiwellPlateTypeEnumNet [GetPlateType](#) ()
Gets the plate type.
- MultiwellPlateTypeEnumNet [GetPlateTypeByHeadstage](#) (uint32_t Headstage)
Gets the plate type.
- void [SetPlateType](#) (MultiwellPlateTypeEnumNet plateType)
Sets the plate type.
- void [SetPlateTypeByHeadstage](#) (uint32_t Headstage, MultiwellPlateTypeEnumNet plateType)
Sets the plate type.
- void [SetPlateMux](#) (uint32_t muxSelection)
Selects a one quarter of the electrodes on a high density Multiwell plate.
- void [SetPlateMuxByHeadstage](#) (uint32_t Headstage, uint32_t muxSelection)
Selects a one quarter of the electrodes on a high density Multiwell plate.
- uint32_t [GetPlateMux](#) ()
Gets the selected quarter of the electrodes on a high density Multiwell plate.
- uint32_t [GetPlateMuxByHeadstage](#) (uint32_t Headstage)
Gets the selected quarter of the electrodes on a high density Multiwell plate.
- bool [IsPlateTypeValid](#) ()
Checks whether the plate type is valid, meaning all pins have contact.
- bool [IsPlateTypeValidByHeadstage](#) (uint32_t Headstage)
Checks whether the plate type is valid, meaning all pins have contact.

Additional Inherited Members

10.70.1 Detailed Description

[CMultiwellDeviceNet](#) is the class to access the HEKA LIH3 device.

10.70.2 Constructor & Destructor Documentation

10.70.2.1 CMultiwellDeviceNet() `CMultiwellDeviceNet ()`

Initializes a new instance of the `CMultiwellDeviceNet` class.

10.70.2.2 ~CMultiwellDeviceNet() `virtual ~CMultiwellDeviceNet () [virtual]`

10.70.2.3 "!CMultiwellDeviceNet() `!CMultiwellDeviceNet ()`

10.70.3 Member Function Documentation

10.70.3.1 ClosePlateClamp() `void ClosePlateClamp ()`

Closes the plate clamp.

10.70.3.2 GetPlateClampLockState() `uint32_t GetPlateClampLockState ()`

Gets the state of the plate clamp lock.

Returns

the state of the plate lock (unlocked/locked)

10.70.3.3 GetPlateClampState() `PlateClampEnumNet GetPlateClampState ()`

Gets the state of the Multiwell plate clamp.

Returns

the state of the plate clamp (open/closed)

10.70.3.4 GetPlateClampStateByHeadstage() `PlateClampEnumNet GetPlateClampStateByHeadstage (uint32_t Headstage)`

Gets the state of the plate

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The plate state

10.70.3.5 GetPlateMux() `uint32_t GetPlateMux ()`

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Returns

the selected quarter

10.70.3.6 GetPlateMuxByHeadstage() `uint32_t GetPlateMuxByHeadstage (
uint32_t Headstage)`

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

the selected quarter

10.70.3.7 GetPlateType() `MultiwellPlateTypeEnumNet GetPlateType ()`

Gets the plate type.

Returns

the plate type

10.70.3.8 GetPlateTypeByHeadstage() `MultiwellPlateTypeEnumNet GetPlateTypeByHeadstage (
uint32_t Headstage)`

Gets the plate type.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

the plate type

10.70.3.9 IsPlateTypeValid() `bool IsPlateTypeValid ()`

Checks whether the plate type is valid, meaning all pins have contact.

Returns

"true" when all pins have contact, otherwise "false".

10.70.3.10 IsPlateTypeValidByHeadstage() `bool IsPlateTypeValidByHeadstage (
uint32_t Headstage)`

Checks whether the plate type is valid, meaning all pins have contact.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

"true" when all pins have contact, otherwise "false".

10.70.3.11 LockPlateClamp() `void LockPlateClamp ()`

Locks the plate clamp.

10.70.3.12 OpenPlateClamp() `void OpenPlateClamp ()`

Opens the plate clamp.

10.70.3.13 SetPlateMux() `void SetPlateMux (
uint32_t muxSelection)`

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>muxSelection</i>	the selected quarter
---------------------	----------------------

10.70.3.14 SetPlateMuxByHeadstage() `void SetPlateMuxByHeadstage (`
 `uint32_t Headstage,`
 `uint32_t muxSelection)`

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>Headstage</i>	The headstage to query.
<i>muxSelection</i>	the selected quarter

10.70.3.15 SetPlateType() `void SetPlateType (`
 `MultiwellPlateTypeEnumNet plateType)`

Sets the plate type.

Parameters

<i>plateType</i>	the plate type
------------------	----------------

10.70.3.16 SetPlateTypeByHeadstage() `void SetPlateTypeByHeadstage (`
 `uint32_t Headstage,`
 `MultiwellPlateTypeEnumNet plateType)`

Sets the plate type.

Parameters

<i>Headstage</i>	The headstage to query.
<i>plateType</i>	the plate type

10.70.3.17 StopPlateClamp() `void StopPlateClamp ()`

Stops the plate clamp movement.


10.70.3.18 UnlockPlateClamp() void UnlockPlateClamp ()

Unlocks the plate clamp.

10.71 CMultiwellOptoStimFunctionNet Class Reference

[CMultiwellOptoStimFunctionNet](#) is the class to access the optical properties of the Multiwell Optostim device

Inheritance diagram for CMultiwellOptoStimFunctionNet:

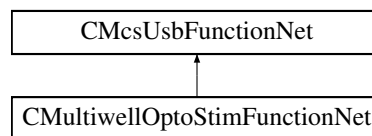**Public Member Functions**

- [CMultiwellOptoStimFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pMultiwellOptoStimFunctionPointerContainer)
Initializes a new instance of the [CMultiwellOptoStimFunctionNet](#) class.
- [CMultiwellOptoStimFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMultiwellOptoStimFunctionNet](#) ()
- [!CMultiwellOptoStimFunctionNet](#) ()
- [uint32_t](#) [GetWaveLengthInNanometer](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetAbsMaxCurrentInMicroAmp](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetMaxDurationHighCurrentInMicroSec](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetMaxDutyCycleHighCurrent](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetPermanentCurrentInMicroAmp](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetColorRgb](#) ([uint16_t](#) channel)
- [String](#)[^] [GetColorStr](#) ([uint16_t](#) channel)
- void [SetWaveLengthInNanometer](#) ([uint16_t](#) channel, [uint32_t](#) WaveLength_nm)
- void [SetAbsMaxCurrentInMicroAmp](#) ([uint16_t](#) channel, [uint32_t](#) AbsoluteMaxCurrent_uA)
- void [SetMaxDurationHighCurrentInMicroSec](#) ([uint16_t](#) channel, [uint32_t](#) AbsoluteMaxDuration_us)
- void [SetMaxDutyCycleHighCurrent](#) ([uint16_t](#) channel, [uint32_t](#) MaxDutyCycleHighCurrent)
- void [SetPermanentCurrentInMicroAmp](#) ([uint16_t](#) channel, [uint32_t](#) PermanentCurrent_uA)
- void [SetColorRgb](#) ([uint16_t](#) channel, [uint32_t](#) ColorRGB)
- void [SetColorStr](#) ([uint16_t](#) channel, [String](#)[^] ColorString)

Additional Inherited Members**10.71.1 Detailed Description**

[CMultiwellOptoStimFunctionNet](#) is the class to access the optical properties of the Multiwell Optostim device

10.71.2 Constructor & Destructor Documentation

10.71.2.1 CMultiwellOptoStimFunctionNet() [1/2] `CMultiwellOptoStimFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pMultiwellOptoStimFunctionPointerContainer)`

Initializes a new instance of the `CMultiwellOptoStimFunctionNet` class.

10.71.2.2 CMultiwellOptoStimFunctionNet() [2/2] `CMultiwellOptoStimFunctionNet (CMcsUsbNet^ mcsusb)`

10.71.2.3 ~CMultiwellOptoStimFunctionNet() `virtual ~CMultiwellOptoStimFunctionNet ()` [virtual]

10.71.2.4 "!CMultiwellOptoStimFunctionNet() `!CMultiwellOptoStimFunctionNet ()`

10.71.3 Member Function Documentation

10.71.3.1 GetAbsMaxCurrentInMicroAmp() `uint32_t GetAbsMaxCurrentInMicroAmp (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

absolute max. current; unit: uA

10.71.3.2 GetColorRgb() `uint32_t GetColorRgb (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

RGB-value of LED color

10.71.3.3 GetColorStr() `String ^ GetColorStr (`
`uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

LED color as string

10.71.3.4 GetMaxDurationHighCurrentInMicroSec() `uint32_t GetMaxDurationHighCurrentInMicroSec (`
`uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. duration the LED can stand the abs. max current; unit: us

10.71.3.5 GetMaxDutyCycleHighCurrent() `uint32_t GetMaxDutyCycleHighCurrent (`
`uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. duty cycle at max. current; unit: 100*%

10.71.3.6 GetPermanentCurrentInMicroAmp() `uint32_t GetPermanentCurrentInMicroAmp (`
`uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. current the LED can stand when always switched on; unit: uA

10.71.3.7 GetWaveLengthInNanometer() `uint32_t GetWaveLengthInNanometer (`
`uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

wavelength of this channel's LEDs; unit: nm

10.71.3.8 SetAbsMaxCurrentInMicroAmp() `void SetAbsMaxCurrentInMicroAmp (`
`uint16_t channel,`
`uint32_t AbsoluteMaxCurrent_uA)`

Parameters

<i>channel</i>	the (analog) channel number
<i>AbsoluteMaxCurrent_uA</i>	absolute max. current; unit: uA

10.71.3.9 SetColorRgb() `void SetColorRgb (`
`uint16_t channel,`
`uint32_t ColorRGB)`

Parameters

<i>channel</i>	the (analog) channel number
<i>ColorRGB</i>	RGB-value of LED color

10.71.3.10 SetColorStr() `void SetColorStr (`
`uint16_t channel,`
`String^ ColorString)`

Parameters

<i>channel</i>	the (analog) channel number
<i>ColorString</i>	LED color as string

10.71.3.11 SetMaxDurationHighCurrentInMicroSec() void SetMaxDurationHighCurrentInMicroSec (
 uint16_t *channel*,
 uint32_t *AbsoluteMaxDuration_us*)

Parameters

<i>channel</i>	the (analog) channel number
<i>AbsoluteMaxDuration_us</i>	max. duration the LED can stand the abs. max current; unit: us

10.71.3.12 SetMaxDutyCycleHighCurrent() void SetMaxDutyCycleHighCurrent (
 uint16_t *channel*,
 uint32_t *MaxDutyCycleHighCurrent*)

Parameters

<i>channel</i>	the (analog) channel number
<i>MaxDutyCycleHighCurrent</i>	max. duty cycle at max. current; unit: 100*%

10.71.3.13 SetPermanentCurrentInMicroAmp() void SetPermanentCurrentInMicroAmp (
 uint16_t *channel*,
 uint32_t *PermanentCurrent_uA*)

Parameters

<i>channel</i>	the (analog) channel number
<i>PermanentCurrent_uA</i>	max. current the LED can stand when always switched on; unit: uA


10.71.3.14 SetWaveLengthInNanometer() void SetWaveLengthInNanometer (
 uint16_t *channel*,
 uint32_t *WaveLength_nm*)

Parameters

<i>channel</i>	the (analog) channel number
<i>WaveLength_nm</i>	wavelength of this channel's LEDs; unit: nm

10.72 CNF_GenDeviceNet Class Reference

Inheritance diagram for CNF_GenDeviceNet:

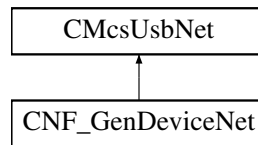

Public Member Functions

- [CNF_GenDeviceNet](#) (void)
- [~CNF_GenDeviceNet](#) (void)
- void [Set_Values](#) (unsigned int frequency, unsigned int amplitude)

Additional Inherited Members

10.72.1 Constructor & Destructor Documentation

10.72.1.1 [CNF_GenDeviceNet\(\)](#) [CNF_GenDeviceNet](#) (
 void)


10.72.1.2 [~CNF_GenDeviceNet\(\)](#) [~CNF_GenDeviceNet](#) (
 void)

10.72.2 Member Function Documentation

10.72.2.1 [Set_Values\(\)](#) void [Set_Values](#) (
 unsigned int *frequency*,
 unsigned int *amplitude*)

10.73 COctoPotDeviceNet Class Reference

Inheritance diagram for COctoPotDeviceNet:

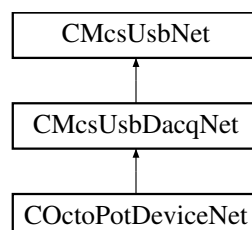

Public Member Functions

- [COctoPotDeviceNet](#) (void)
- [COctoPotDeviceNet](#) ([OnChannelData](#)^ channelData, [OnError](#)^ error)
- uint32_t [SetOutputRate](#) (uint32_t rate)
- uint32_t [SetBathclamp](#) (unsigned int block, bool enable)
- uint32_t [SetDacValue](#) (int channel, int value)
- uint32_t [SetDacAutoControl](#) (unsigned int channel)
- uint32_t [SetPidParameter](#) (unsigned int channel, int const_p, int const_i, int shift_p, int shift_i)
- uint32_t [SetRampParameter](#) (unsigned int channel, int start, int min, int max, int slope, int slope2, int pause, unsigned int samples)
- uint32_t [RampStart](#) (int channelmap)
- uint32_t [SetSineParameter](#) (unsigned int channel, int amplitude)
- uint32_t [SineStart](#) (int channelmap)
- uint32_t [SetPatternListEntry](#) (unsigned int channel, unsigned int position, unsigned int duration, int value)
- uint32_t [PatternListStart](#) (int channelmap)
- uint32_t [SetAdcOffset](#) (unsigned int channel, int offset)
- uint32_t [SetDacOffset](#) (unsigned int channel, int offset)
- uint32_t [ResetAdcOffset](#) (unsigned int channel)
- uint32_t [ResetDacOffset](#) (unsigned int channel)
- uint32_t [BurnAdcOffset](#) ()
- uint32_t [BurnDacOffset](#) ()
- uint32_t [GetAdcOffset](#) (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- uint32_t [GetDacOffset](#) (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- uint32_t [SetAmplificationSwitch](#) (unsigned int channel, unsigned int state)
- uint32_t [SetChannelSwitch](#) (unsigned int channel, unsigned int state)
- uint32_t [SetNumberOfChannels](#) (unsigned int NumberOfChannels)
- uint32_t [EnableDigitalIn](#) (bool enable)
- uint32_t [EnableTimestamp](#) (bool enable)
- uint32_t [EnableChecksum](#) (bool enable)

Additional Inherited Members**10.73.1 Constructor & Destructor Documentation**

10.73.1.1 COctoPotDeviceNet() [1/2] [COctoPotDeviceNet](#) (void)

10.73.1.2 COctoPotDeviceNet() [2/2] [COctoPotDeviceNet](#) ([OnChannelData](#)^ channelData, [OnError](#)^ error)

10.73.2 Member Function Documentation

10.73.2.1 BurnAdcOffset() uint32_t BurnAdcOffset ()

10.73.2.2 BurnDacOffset() uint32_t BurnDacOffset ()

10.73.2.3 EnableChecksum() uint32_t EnableChecksum (
bool *enable*)

10.73.2.4 EnableDigitalIn() uint32_t EnableDigitalIn (
bool *enable*)

10.73.2.5 EnableTimestamp() uint32_t EnableTimestamp (
bool *enable*)

10.73.2.6 GetAdcOffset() uint32_t GetAdcOffset (
unsigned int *channel*,
[System::Runtime::InteropServices::Out] int ^ *offset*)

10.73.2.7 GetDacOffset() uint32_t GetDacOffset (
unsigned int *channel*,
[System::Runtime::InteropServices::Out] int ^ *offset*)

10.73.2.8 PatternListStart() uint32_t PatternListStart (
int *channelmap*)

10.73.2.9 RampStart() uint32_t RampStart (
int *channelmap*)

10.73.2.10 ResetAdcOffset() uint32_t ResetAdcOffset (
unsigned int *channel*)

10.73.2.11 ResetDacOffset() uint32_t ResetDacOffset (
 unsigned int *channel*)

10.73.2.12 SetAdcOffset() uint32_t SetAdcOffset (
 unsigned int *channel*,
 int *offset*)

10.73.2.13 SetAmplificationSwitch() uint32_t SetAmplificationSwitch (
 unsigned int *channel*,
 unsigned int *state*)

10.73.2.14 SetBathclamp() uint32_t SetBathclamp (
 unsigned int *block*,
 bool *enable*)

10.73.2.15 SetChannelSwitch() uint32_t SetChannelSwitch (
 unsigned int *channel*,
 unsigned int *state*)

10.73.2.16 SetDacAutoControl() uint32_t SetDacAutoControl (
 unsigned int *channel*)

10.73.2.17 SetDacOffset() uint32_t SetDacOffset (
 unsigned int *channel*,
 int *offset*)

10.73.2.18 SetDacValue() uint32_t SetDacValue (
 int *channel*,
 int *value*)

10.73.2.19 SetNumberOfChannels() uint32_t SetNumberOfChannels (
 unsigned int *NumberOfChannels*)

10.73.2.20 SetOutputRate() uint32_t SetOutputRate (
 uint32_t *rate*)

10.73.2.21 SetPatternListEntry() uint32_t SetPatternListEntry (
 unsigned int *channel*,
 unsigned int *position*,
 unsigned int *duration*,
 int *value*)

10.73.2.22 SetPidParameter() uint32_t SetPidParameter (
 unsigned int *channel*,
 int *const_p*,
 int *const_i*,
 int *shift_p*,
 int *shift_i*)


10.73.2.23 SetRampParameter() uint32_t SetRampParameter (
 unsigned int *channel*,
 int *start*,
 int *min*,
 int *max*,
 int *slope*,
 int *slope2*,
 int *pause*,
 unsigned int *samples*)

10.73.2.24 SetSineParameter() uint32_t SetSineParameter (
 unsigned int *channel*,
 int *amplitude*)

10.73.2.25 SineStart() uint32_t SineStart (
 int *channelmap*)

10.74 COkuvisionStimulatorDeviceNet Class Reference

Inheritance diagram for COkuvisionStimulatorDeviceNet:

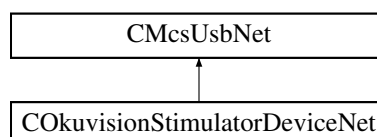

Public Member Functions

- [COkuvisionStimulatorDeviceNet](#) (void)
- [~COkuvisionStimulatorDeviceNet](#) (void)
- void [SetPulseform](#) (int channel, int current, int pulsewidth, int periode, int duration)
- void [GetPulseform](#) (int channel, [System::Runtime::InteropServices::Out] int% current, [System::Runtime::InteropServices::Out] int% pulsewidth, [System::Runtime::InteropServices::Out] int% periode, [System::Runtime::InteropServices::Out] int% duration)
- void [SetMaxPower](#) (int channel, int power)
- int [GetMaxPower](#) (int channel)
- void [SetMaxVoltage](#) (int channel, int voltage)
- int [GetMaxVoltage](#) (int channel)
- void [SetCheckVoltage](#) (int channel, int voltage)
- int [GetCheckVoltage](#) (int channel)
- int [GetVoltage](#) (int channel)
- void [SetDACOffset](#) (int channel, int part, int offset)
- int [GetDACOffset](#) (int channel, int part)
- void [SetRTC](#) (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)
- void [GetRTC](#) ([System::Runtime::InteropServices::Out] uint8_t% year, [System::Runtime::InteropServices::Out] uint8_t% month, [System::Runtime::InteropServices::Out] uint8_t% day, [System::Runtime::InteropServices::Out] uint8_t% hour, [System::Runtime::InteropServices::Out] uint8_t% minute, [System::Runtime::InteropServices::Out] uint8_t% second)
- void [SetRTC](#) (DateTime timestamp)
- DateTime [GetRTC](#) ()
- void [GetStimulatorStatus](#) ([System::Runtime::InteropServices::Out] int% startstop, [System::Runtime::InteropServices::Out] int% last_error, [System::Runtime::InteropServices::Out] int% battery_status)
- void [SetCurrentFactor](#) (int channel, int factor)
- int [GetCurrentFactor](#) (int channel)

Additional Inherited Members**10.74.1 Constructor & Destructor Documentation**

10.74.1.1 COkuvisionStimulatorDeviceNet() [COkuvisionStimulatorDeviceNet](#) (
 void)

10.74.1.2 ~COkuvisionStimulatorDeviceNet() [~COkuvisionStimulatorDeviceNet](#) (
 void)

10.74.2 Member Function Documentation

10.74.2.1 GetCheckVoltage() int GetCheckVoltage (
 int channel)

10.74.2.2 GetCurrentFactor() `int GetCurrentFactor (`
 `int channel)`

10.74.2.3 GetDACOffset() `int GetDACOffset (`
 `int channel,`
 `int part)`

10.74.2.4 GetMaxPower() `int GetMaxPower (`
 `int channel)`

10.74.2.5 GetMaxVoltage() `int GetMaxVoltage (`
 `int channel)`

10.74.2.6 GetPulseform() `void GetPulseform (`
 `int channel,`
 `[System::Runtime::InteropServices::Out] int% current,`
 `[System::Runtime::InteropServices::Out] int% pulsewidth,`
 `[System::Runtime::InteropServices::Out] int% periode,`
 `[System::Runtime::InteropServices::Out] int% duration)`

10.74.2.7 GetRTC() [1/2] `DateTime GetRTC ()`

10.74.2.8 GetRTC() [2/2] `void GetRTC (`
 `[System::Runtime::InteropServices::Out] uint8_t% year,`
 `[System::Runtime::InteropServices::Out] uint8_t% month,`
 `[System::Runtime::InteropServices::Out] uint8_t% day,`
 `[System::Runtime::InteropServices::Out] uint8_t% hour,`
 `[System::Runtime::InteropServices::Out] uint8_t% minute,`
 `[System::Runtime::InteropServices::Out] uint8_t% second)`

10.74.2.9 GetStimulatorStatus() `void GetStimulatorStatus (`
 `[System::Runtime::InteropServices::Out] int% startstop,`
 `[System::Runtime::InteropServices::Out] int% last_error,`
 `[System::Runtime::InteropServices::Out] int% battery_status)`

10.74.2.10 GetVoltage() `int GetVoltage (`
`int channel)`

10.74.2.11 SetCheckVoltage() `void SetCheckVoltage (`
`int channel,`
`int voltage)`

10.74.2.12 SetCurrentFactor() `void SetCurrentFactor (`
`int channel,`
`int factor)`

10.74.2.13 SetDACOffset() `void SetDACOffset (`
`int channel,`
`int part,`
`int offset)`

10.74.2.14 SetMaxPower() `void SetMaxPower (`
`int channel,`
`int power)`

10.74.2.15 SetMaxVoltage() `void SetMaxVoltage (`
`int channel,`
`int voltage)`

10.74.2.16 SetPulseform() `void SetPulseform (`
`int channel,`
`int current,`
`int pulsewidth,`
`int periode,`
`int duration)`


10.74.2.17 SetRTC() `[1/2] void SetRTC (`
`DateTime timestamp)`

10.74.2.18 SetRTC() [2/2] `void SetRTC (`
`uint8_t year,`
`uint8_t month,`
`uint8_t day,`
`uint8_t hour,`
`uint8_t minute,`
`uint8_t second)`

10.75 CPatchServerDeviceNet Class Reference

[CPatchServerDeviceNet](#) is the class to control the MCS PatchServer device

Inheritance diagram for CPatchServerDeviceNet:

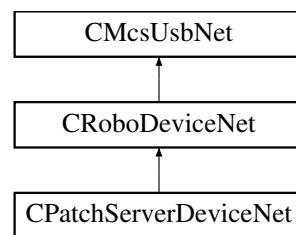

Public Member Functions

- [CPatchServerDeviceNet](#) (void)

Properties

- [CMcsBus_SensorNet](#)[^] [Sensor](#) [get]

Additional Inherited Members

10.75.1 Detailed Description

[CPatchServerDeviceNet](#) is the class to control the MCS PatchServer device

10.75.2 Constructor & Destructor Documentation


10.75.2.1 CPatchServerDeviceNet() [CPatchServerDeviceNet](#) (
`void`)

10.75.3 Property Documentation

10.75.3.1 Sensor `CMcsBus_SensorNet`^ Sensor [get]

10.76 CPathIdentDeviceNet Class Reference

Inheritance diagram for CPathIdentDeviceNet:

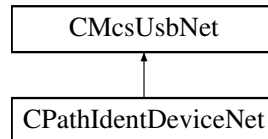

Public Member Functions

- `CPathIdentDeviceNet` (void)
- `~CPathIdentDeviceNet` (void)
- void `Set_Values` (unsigned int frequency, unsigned int amplitude)
- void `Measure` ([System::Runtime::InteropServices::Out] unsigned int% phase, [System::Runtime::InteropServices::Out] unsigned int% amplitude)

Additional Inherited Members

10.76.1 Constructor & Destructor Documentation

10.76.1.1 CPathIdentDeviceNet() `CPathIdentDeviceNet` (void)

10.76.1.2 ~CPathIdentDeviceNet() `~CPathIdentDeviceNet` (void)


10.76.2 Member Function Documentation

10.76.2.1 Measure() void Measure ([System::Runtime::InteropServices::Out] unsigned int% *phase*, [System::Runtime::InteropServices::Out] unsigned int% *amplitude*)

10.76.2.2 Set_Values() `void Set_Values (`
 `unsigned int frequency,`
 `unsigned int amplitude)`

10.77 CPedoterDeviceNet Class Reference

Inheritance diagram for CPedoterDeviceNet:

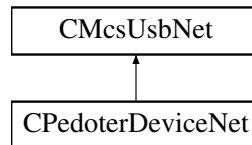

Public Member Functions

- [CPedoterDeviceNet](#) ()
Initializes a new instance of the [CPedoterDeviceNet](#) class.
- virtual [~CPedoterDeviceNet](#) ()
- [!CPedoterDeviceNet](#) ()
- uint32_t [GetCommand](#) (uint16_t Argument)
Get value from the pedoter device
- void [SetCommand](#) (uint16_t Argument, uint32_t pData)
Set value on the pedoter device

Additional Inherited Members

10.77.1 Detailed Description

10.77.2 Constructor & Destructor Documentation

10.77.2.1 CPedoterDeviceNet() [CPedoterDeviceNet](#) ()

Initializes a new instance of the [CPedoterDeviceNet](#) class.

10.77.2.2 ~CPedoterDeviceNet() `virtual ~CPedoterDeviceNet ()` [virtual]

10.77.2.3 "!CPedoterDeviceNet() [!CPedoterDeviceNet](#) ()

10.77.3 Member Function Documentation

10.77.3.1 GetCommand() `uint32_t GetCommand (`
 `uint16_t Argument)`

Get value from the pedoter device

Parameters

<i>Argument</i>	argument
-----------------	----------

Returns

value

10.77.3.2 SetCommand() void SetCommand (
 uint16_t *Argument*,
 uint32_t *pData*)

Set value on the pedoter device


Parameters

<i>Argument</i>	argument
<i>pData</i>	value

10.78 CPeristalticPumpDeviceNet Class Reference

[CPeristalticPumpDeviceNet](#) is the class to control a Persistaltic Pump.

Inheritance diagram for CPeristalticPumpDeviceNet:

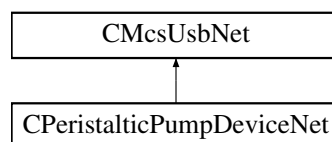

Public Member Functions

- [CPeristalticPumpDeviceNet](#) (void)
 Initialize a new instance of the [CPeristalticPumpDeviceNet](#) class.
- [~CPeristalticPumpDeviceNet](#) (void)

Properties

- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]

Additional Inherited Members

10.78.1 Detailed Description

[CPeristalticPumpDeviceNet](#) is the class to control a Persistaltic Pump.

10.78.2 Constructor & Destructor Documentation

10.78.2.1 CPeristalticPumpDeviceNet() `CPeristalticPumpDeviceNet (void)`

Initialize a new instance of the `CPeristalticPumpDeviceNet` class.


10.78.2.2 ~CPeristalticPumpDeviceNet() `~CPeristalticPumpDeviceNet (void)`

10.78.3 Property Documentation

10.78.3.1 McsBus_MotorControl `CMcsBus_MotorControlNet^ McsBus_MotorControl [get]`

10.79 CPgaDeviceNet Class Reference

Inheritance diagram for CPgaDeviceNet:

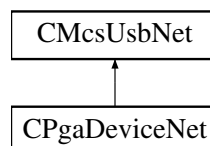

Public Member Functions

- `CPgaDeviceNet ()`
- `~CPgaDeviceNet ()`
- `uint32_t GetNumFrequencyRanges ([System::Runtime::InteropServices::Out]int% numRanges)`
- `uint32_t GetFrequencyRange (int rangeIndex, [System::Runtime::InteropServices::Out]int% low, [System::Runtime::InteropServices::Out]int% high, [System::Runtime::InteropServices::Out]int% channels, [System::Runtime::InteropServices::Out]int% gain)`
- `uint32_t GetNumAmplifications ([System::Runtime::InteropServices::Out]int% number)`
- `uint32_t GetAmplification (int index, [System::Runtime::InteropServices::Out]int% amplification, [System::Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)`
- `uint32_t DefineNumFrequencyRanges (int rnum)`
- `uint32_t DefineFrequencyRange (int index, int low, int high, int channels, int gain)`
- `uint32_t DefineNumAmplifications (int number)`
- `uint32_t DefineAmplification (int index, int amplification, int poti1, int poti2)`
- `uint32_t SetGain (int channel, int Gain, int poti1, int poti2)`
- `uint32_t GetGain (int channel, [System::Runtime::InteropServices::Out]int% Gain, [System::Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)`
- `uint32_t ApplyGains ()`

Additional Inherited Members

10.79.1 Constructor & Destructor Documentation

10.79.1.1 CPgaDeviceNet() `CPgaDeviceNet ()`

10.79.1.2 ~CPgaDeviceNet() `~CPgaDeviceNet ()`

10.79.2 Member Function Documentation

10.79.2.1 ApplyGains() `uint32_t ApplyGains ()`

10.79.2.2 DefineAmplification() `uint32_t DefineAmplification (`
 `int index,`
 `int amplification,`
 `int poti1,`
 `int poti2)`

10.79.2.3 DefineFrequencyRange() `uint32_t DefineFrequencyRange (`
 `int index,`
 `int low,`
 `int high,`
 `int channels,`
 `int gain)`

10.79.2.4 DefineNumAmplifications() `uint32_t DefineNumAmplifications (`
 `int number)`

10.79.2.5 DefineNumFrequencyRanges() `uint32_t DefineNumFrequencyRanges (`
 `int rnum)`

10.79.2.6 GetAmplification() uint32_t GetAmplification (
 int *index*,
 [System::Runtime::InteropServices::Out] int% *amplification*,
 [System::Runtime::InteropServices::Out] int% *potil*,
 [System::Runtime::InteropServices::Out] int% *poti2*)

10.79.2.7 GetFrequencyRange() uint32_t GetFrequencyRange (
 int *rangeIndex*,
 [System::Runtime::InteropServices::Out] int% *low*,
 [System::Runtime::InteropServices::Out] int% *high*,
 [System::Runtime::InteropServices::Out] int% *channels*,
 [System::Runtime::InteropServices::Out] int% *gain*)

10.79.2.8 GetGain() uint32_t GetGain (
 int *channel*,
 [System::Runtime::InteropServices::Out] int% *Gain*,
 [System::Runtime::InteropServices::Out] int% *potil*,
 [System::Runtime::InteropServices::Out] int% *poti2*)

10.79.2.9 GetNumAmplifications() uint32_t GetNumAmplifications (
 [System::Runtime::InteropServices::Out] int% *number*)


10.79.2.10 GetNumFrequencyRanges() uint32_t GetNumFrequencyRanges (
 [System::Runtime::InteropServices::Out] int% *numRanges*)

10.79.2.11 SetGain() uint32_t SetGain (
 int *channel*,
 int *Gain*,
 int *potil*,
 int *poti2*)

10.80 CPositionIIDeviceNet Class Reference

[CPositionIIDeviceNet](#) is the class to control PositionII devices

Inheritance diagram for CPositionIIDeviceNet:

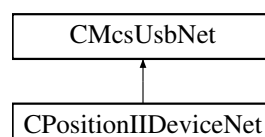

Public Member Functions

- [CPositionIIDeviceNet](#) ()
Initializes a new instance of the [CPositionIIDeviceNet](#) class.
- virtual [~CPositionIIDeviceNet](#) ()
- [!CPositionIIDeviceNet](#) ()
- uint32_t [GetCoilCommunication](#) (uint16_t coil)
get if the communication to the coil is working
- uint32_t [GetOnOff](#) (uint16_t coil)
get if the coil is switched on/off
- void [SwitchOnOff](#) (uint16_t coil, uint32_t on)
switched the coil on of
- uint32_t [GetImplantState](#) (uint16_t coil)
gets the implantat state
- uint32_t [GetImplantCurrentSetpoint](#) (uint16_t coil)
sets the implant current setpoint
- void [SetImplantCurrentSetpoint](#) (uint16_t coil, uint32_t current)
gets the implant current setpoint

Properties

- [CRFFunctionNet](#)[^] [RFFunction](#) [get]

Additional Inherited Members**10.80.1 Detailed Description**

[CPositionIIDeviceNet](#) is the class to control PositionII devices

10.80.2 Constructor & Destructor Documentation**10.80.2.1 CPositionIIDeviceNet()** [CPositionIIDeviceNet](#) ()

Initializes a new instance of the [CPositionIIDeviceNet](#) class.

10.80.2.2 ~CPositionIIDeviceNet() virtual [~CPositionIIDeviceNet](#) () [virtual]**10.80.2.3 !CPositionIIDeviceNet()** [!CPositionIIDeviceNet](#) ()**10.80.3 Member Function Documentation****10.80.3.1 GetCoilCommunication()** uint32_t [GetCoilCommunication](#) (uint16_t coil)

get if the communication to the coil is working

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

is communicating

10.80.3.2 GetImplantCurrentSetpoint() uint32_t GetImplantCurrentSetpoint (
 uint16_t *coil*)

sets the implant current setpoint

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the current

10.80.3.3 GetImplantState() uint32_t GetImplantState (
 uint16_t *coil*)

gets the implantat state

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

10.80.3.4 GetOnOff() uint32_t GetOnOff (
 uint16_t *coil*)

get if the coil is switched on/off

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

0 = off, 1 = on

10.80.3.5 SetImplantCurrentSetpoint() `void SetImplantCurrentSetpoint (`
 `uint16_t coil,`
 `uint32_t current)`

gets the implant current setpoint

Parameters

<i>coil</i>	the coil
<i>current</i>	the current

10.80.3.6 SwitchOnOff() `void SwitchOnOff (`
 `uint16_t coil,`
 `uint32_t on)`

switched the coild on of

Parameters

<i>coil</i>	the coil
<i>on</i>	0 = off, 1 = on


10.80.4 Property Documentation

10.80.4.1 RFFunction `CRFFunctionNet^ RFFunction [get]`

10.81 CPositionImpDeviceNet Class Reference

[CPositionImpDeviceNet](#) is the class to access the Position/Imp devices

Inheritance diagram for CPositionImpDeviceNet:

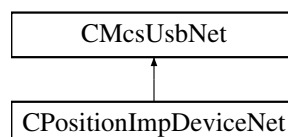

Public Member Functions

- [CPositionImpDeviceNet](#) ()
Initializes a new instance of the [CPositionImpDeviceNet](#) class.
- virtual [~CPositionImpDeviceNet](#) ()
- [!CPositionImpDeviceNet](#) ()
- void [ConnectImp](#) (uint32_t id)
Connect to a Imp device with a certain ID
- uint32_t [ConnectedImp](#) ()
The ID of the connected Imp device
- int32_t [GetRFFrequency](#) ()
Gets currently used RF frequency
- void [SetRFFrequency](#) (int32_t frequency)
Sets the current RF frequency
- uint32_t [GetDeviceList](#) (int32_t index)
Gets the device list
- void [SetDeviceList](#) (int32_t index, uint32_t id)
Sets the device list
- uint32_t [GetImpId](#) ()
Gets the ID of the impedance measure device
- void [SetImpId](#) (uint32_t id)
Sets the ID of the impedance measure device

Additional Inherited Members**10.81.1 Detailed Description**

[CPositionImpDeviceNet](#) is the class to access the Position/Imp devices

10.81.2 Constructor & Destructor Documentation**10.81.2.1 CPositionImpDeviceNet()** [CPositionImpDeviceNet](#) ()

Initializes a new instance of the [CPositionImpDeviceNet](#) class.

10.81.2.2 ~CPositionImpDeviceNet() virtual [~CPositionImpDeviceNet](#) () [virtual]**10.81.2.3 "!CPositionImpDeviceNet()** [!CPositionImpDeviceNet](#) ()

10.81.3 Member Function Documentation

10.81.3.1 ConnectedImp() `uint32_t ConnectedImp ()`

The ID of the connected Imp device

Returns

The ID

10.81.3.2 ConnectImp() `void ConnectImp (uint32_t id)`

Connect to a Imp device with a certain ID

Parameters

<i>id</i>	The ID
-----------	--------

10.81.3.3 GetDeviceList() `uint32_t GetDeviceList (int32_t index)`

Gets the device list

Parameters

<i>index</i>	the index
--------------	-----------

Returns

the ID

10.81.3.4 GetImpId() `uint32_t GetImpId ()`

Gets the ID of the impedance measure device

Returns

the ID

10.81.3.5 GetRFFrequency() `int32_t GetRFFrequency ()`

Gets currently used RF frequency

Returns

The frequency

10.81.3.6 SetDeviceList() `void SetDeviceList (
int32_t index,
uint32_t id)`

Sets the device list

Parameters

<i>index</i>	the index
<i>id</i>	the ID

10.81.3.7 SetImpId() `void SetImpId (
uint32_t id)`

Sets the ID of the impedance measure device

Parameters

<i>id</i>	the ID
-----------	--------

10.81.3.8 SetRFFrequency() `void SetRFFrequency (
int32_t frequency)`


Sets the current RF frequency

Parameters

<i>frequency</i>	The frequency
------------------	---------------

10.82 CPPCDeviceNet Class Reference

Inheritance diagram for CPPCDeviceNet:

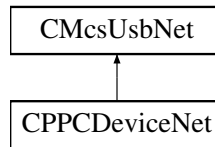

Public Member Functions

- [CPPCDeviceNet](#) (void)

Properties

- [CPPCFunctionNet](#)[^] [PPCFunction](#) [get]
- [CMcsBusNet](#)[^] [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]
- [CMcsBus_SensorNet](#)[^] [McsBus_Sensor](#) [get]

Additional Inherited Members

10.82.1 Constructor & Destructor Documentation

10.82.1.1 [CPPCDeviceNet](#)() [CPPCDeviceNet](#) (
void)

10.82.2 Property Documentation

10.82.2.1 [McsBus](#) [CMcsBusNet](#)[^] [McsBus](#) [get]

10.82.2.2 [McsBus_MotorControl](#) [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]


10.82.2.3 [McsBus_Sensor](#) [CMcsBus_SensorNet](#)[^] [McsBus_Sensor](#) [get]

10.82.2.4 [PPCFunction](#) [CPPCFunctionNet](#)[^] [PPCFunction](#) [get]

10.83 CPPCFunctionNet Class Reference

[CPPCFunctionNet](#) is the class to access the PPC (high precision Patch Peristaltic patch Pump

Inheritance diagram for CPPCFunctionNet:

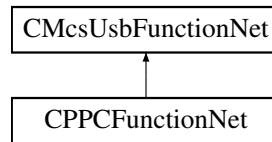

Public Member Functions

- [CPPCFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pPPCFunctionPointer↔ Container)
Initializes a new instance of the [CPPCFunctionNet](#) class.
- [CPPCFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CPPCFunctionNet](#) ()
- [!CPPCFunctionNet](#) ()
- int [GetPumpSpeedUnit](#) (uint16_t channel)
Reads the Pump Speed Unit
- void [SetPumpSpeedUnit](#) (uint16_t channel, int SpeedUnit)
Writes the Pump Speed Unit
- PP_Pump_Mode_Type_EnumNet [GetPumpModeType](#) (uint16_t channel)
Reads the Pump Mode Type.
- void [SetPumpModeType](#) (uint16_t channel, PP_Pump_Mode_Type_EnumNet PumpMode)
Writes the config string from the device.
- void [GetAnalogVoltageRange](#) (uint16_t channel, [System::Runtime::InteropServices::Out]uint16_t% min_↔ voltage, [System::Runtime::InteropServices::Out]uint16_t% max_voltage)
Reads the Analog Input Voltage Range
- void [SetAnalogVoltageRange](#) (uint16_t channel, uint16_t min_voltage, uint16_t max_voltage)
Writes the Analog Input Voltage Range
- void [GetPressureRange](#) (uint16_t channel, [System::Runtime::InteropServices::Out]int32_t% lower_↔ pressure, [System::Runtime::InteropServices::Out]int32_t% upper_pressure)
Get the pressure range that is used between the analog voltage or the digital states
- void [SetPressureRange](#) (uint16_t channel, int32_t lower_pressure, int32_t upper_pressure)
Get the pressure range that is used between the analog voltage or the digital states
- uint16_t [GetSupplyVoltage](#) ()
Reads the current supply voltage in mV
- uint16_t [GetAnalogVoltage](#) (uint16_t channel)
Reads the current analog voltage
- uint16_t [GetDigitalIn](#) (uint16_t channel)
Reads the digital input state
- int [GetValveActive](#) (uint16_t valve)
Gets the valve active/inactive state
- void [SetValveActive](#) (uint16_t valve, int valveActive)
Sets the valve active/inactive state
- void [SetPressureOffset](#) ()
Sets the pressure offset
- void [LoadPressure](#) (int32_t pressure, uint32_t options)

- Loads the reservoir with a pressure*
- void [IsBusy](#) ([System::Runtime::InteropServices::Out]int16_t% task, [System::Runtime::InteropServices::Out]int16_t% wait)
- Is the PPC busy with a task*
- void [FirePressurePulse](#) (int32_t duration, int32_t nextpressure)
- Fire a pressure pulse from the reservoir*
- int32_t [MeasureReservoir](#) ()
- Measures the reservoir pressure*

Additional Inherited Members

10.83.1 Detailed Description

[CPPCFunctionNet](#) is the class to access the PPC (high precision Patch Peristaltic patch Pump

10.83.2 Constructor & Destructor Documentation

10.83.2.1 CPPCFunctionNet() [1/2] [CPPCFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pPPCFunctionPointerContainer)

Initializes a new instance of the [CPPCFunctionNet](#) class.

10.83.2.2 CPPCFunctionNet() [2/2] [CPPCFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.83.2.3 ~CPPCFunctionNet() virtual ~[CPPCFunctionNet](#) () [virtual]

10.83.2.4 "!CPPCFunctionNet() ![CPPCFunctionNet](#) ()

10.83.3 Member Function Documentation

10.83.3.1 FirePressurePulse() void [FirePressurePulse](#) (
int32_t duration,
int32_t nextpressure)

Fire a pressure pulse from the reservoir

Parameters

<i>duration</i>	The pulse duration (valves open)
<i>nextpressure</i>	The next pressure

10.83.3.2 GetAnalogVoltage() `uint16_t GetAnalogVoltage (`
`uint16_t channel)`

Reads the current analog voltage

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Analog Voltage

10.83.3.3 GetAnalogVoltageRange() `void GetAnalogVoltageRange (`
`uint16_t channel,`
`[System::Runtime::InteropServices::Out] uint16_t% min_voltage,`
`[System::Runtime::InteropServices::Out] uint16_t% max_voltage)`

Reads the Analog Input Voltage Range

Parameters

<i>channel</i>	The Channel Number
<i>min_voltage</i>	The voltage that should be seen as the minimum voltage
<i>max_voltage</i>	The voltage that should be seen as the maximum voltage

10.83.3.4 GetDigitalIn() `uint16_t GetDigitalIn (`
`uint16_t channel)`

Reads the digital input state

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Digital State

10.83.3.5 GetPressureRange() void GetPressureRange (
 uint16_t *channel*,
 [System::Runtime::InteropServices::Out] int32_t% *lower_pressure*,
 [System::Runtime::InteropServices::Out] int32_t% *upper_pressure*)

Get the pressure range that is used between the analog voltage or the digital states

Parameters

<i>channel</i>	The Channel Number
<i>lower_pressure</i>	The lower border of the pressure range
<i>upper_pressure</i>	The upper border of the pressure range

10.83.3.6 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType (
 uint16_t *channel*)

Reads the Pump Mode Type.

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Pump Mode Type.

10.83.3.7 GetPumpSpeedUnit() int GetPumpSpeedUnit (
 uint16_t *channel*)

Reads the Pump Speed Unit

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Speed Unit

10.83.3.8 GetSupplyVoltage() `uint16_t GetSupplyVoltage ()`

Reads the current supply voltage in mV

Returns

The supply voltage

10.83.3.9 GetValveActive() `int GetValveActive (
 uint16_t valve)`

Gets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve state

10.83.3.10 IsBusy() `void IsBusy (
 [System::Runtime::InteropServices::Out] int16_t% task,
 [System::Runtime::InteropServices::Out] int16_t% wait)`

Is the PPC busy with a task

Parameters

<i>task</i>	The task state
<i>wait</i>	The wait state

10.83.3.11 LoadPressure() `void LoadPressure (
 int32_t pressure,
 uint32_t options)`

Loads the reservoir with a pressure

Parameters

<i>pressure</i>	The pressure
<i>options</i>	The options: end with 0=regulate on patch 1=regulate on reservoir

10.83.3.12 MeasureReservoir() `int32_t MeasureReservoir ()`

Measures the reservoir pressure

Returns

The pressure

10.83.3.13 SetAnalogVoltageRange() `void SetAnalogVoltageRange (`
`uint16_t channel,`
`uint16_t min_voltage,`
`uint16_t max_voltage)`

Writes the Analog Input Voltage Range

Parameters

<i>channel</i>	The Channel Number
<i>min_voltage</i>	The voltage that should be seen as the minimum voltage
<i>max_voltage</i>	The voltage that should be seen as the maximum voltage

10.83.3.14 SetPressureOffset() `void SetPressureOffset ()`

Sets the pressure offset

10.83.3.15 SetPressureRange() `void SetPressureRange (`
`uint16_t channel,`
`int32_t lower_pressure,`
`int32_t upper_pressure)`

Get the pressure range that is used between the analog voltage or the digital states

Parameters

<i>channel</i>	The Channel Number
<i>lower_pressure</i>	The lower border of the pressure range
<i>upper_pressure</i>	The upper border of the pressure range

10.83.3.16 SetPumpModeType() `void SetPumpModeType (`
 `uint16_t channel,`
 `PP_Pump_Mode_Type_EnumNet PumpMode)`

Writes the config string from the device.

Parameters

<i>channel</i>	The Channel Number
<i>PumpMode</i>	The Pump Mode Type.

10.83.3.17 SetPumpSpeedUnit() `void SetPumpSpeedUnit (`
 `uint16_t channel,`
 `int SpeedUnit)`

Writes the Pump Speed Unit

Parameters

<i>channel</i>	The Channel Number
<i>SpeedUnit</i>	The Speed Unit

10.83.3.18 SetValveActive() `void SetValveActive (`
 `uint16_t valve,`
 `int valveActive)`


Sets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
<i>valveActive</i>	The valve state

10.84 CPPS_DeviceNet Class Reference

Inheritance diagram for CPPS_DeviceNet:

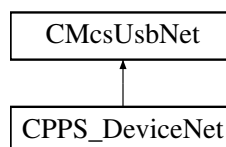

Public Member Functions

- [CPPS_DeviceNet](#) (void)

Properties

- [CPPS_FunctionNet](#)^ [PPS_Function](#) [get]
- [CMcsBusNet](#)^ [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)^ [McsBus_MotorControl](#) [get]
- [CMcsBus_SensorNet](#)^ [McsBus_Sensor](#) [get]

Additional Inherited Members

10.84.1 Constructor & Destructor Documentation

10.84.1.1 [CPPS_DeviceNet\(\)](#) [CPPS_DeviceNet](#) (
void)

10.84.2 Property Documentation

10.84.2.1 [McsBus](#) [CMcsBusNet](#)^ [McsBus](#) [get]


10.84.2.2 [McsBus_MotorControl](#) [CMcsBus_MotorControlNet](#)^ [McsBus_MotorControl](#) [get]

10.84.2.3 [McsBus_Sensor](#) [CMcsBus_SensorNet](#)^ [McsBus_Sensor](#) [get]

10.84.2.4 [PPS_Function](#) [CPPS_FunctionNet](#)^ [PPS_Function](#) [get]

10.85 CPPS_FunctionNet Class Reference

Inheritance diagram for CPPS_FunctionNet:

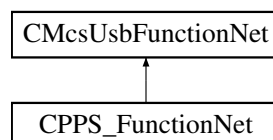

Public Member Functions

- [CPPS_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] cPPS_FunctionPointerContainer)
- [CPPS_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetPumpMaxSpeed](#) (unsigned short index, unsigned short maxspeed)
- unsigned short [GetPumpMaxSpeed](#) (unsigned short index)
- void [SetPumpSpeedUnit](#) (unsigned short index, int speedunit)
- int [GetPumpSpeedUnit](#) (unsigned short index)
- void [SetPumpModeType](#) (unsigned short index, PP_Pump_Mode_Type_EnumNet type)
- PP_Pump_Mode_Type_EnumNet [GetPumpModeType](#) (unsigned short index)
- void [SetPumpCouple](#) (unsigned int i)
- unsigned int [GetPumpCouple](#) ()
- void [SetPumpEnableSpeedRatio](#) (unsigned int enable)
- unsigned int [GetPumpEnableSpeedRatio](#) ()
- void [SetPumpManualOnOff](#) (unsigned short index, unsigned int onoff)
- unsigned int [GetPumpManualOnOff](#) (unsigned short index)
- void [SetPumpFunctionSpeeds](#) (unsigned short index, short offspeed, short onspeed)
- void [GetPumpFunctionSpeeds](#) (unsigned short index, [System::Runtime::InteropServices::Out]short% offspeed, [System::Runtime::InteropServices::Out]short% onspeed)
- void [SetPumpSpeedRatio](#) (int ratio)
- int [GetPumpSpeedRatio](#) ()
- void [SetPumpFastOnOff](#) (unsigned short index, unsigned int onoff)
- unsigned int [GetPumpFastOnOff](#) (unsigned short index)
- void [SetPumpFastSpeed](#) (unsigned short index, short fastspeed)
- short [GetPumpFastSpeed](#) (unsigned short index)
- void [SetAnalogVoltages](#) (unsigned short index, unsigned short minvoltage, unsigned short maxvoltage)
- void [GetAnalogVoltages](#) (unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% minvoltage, [System::Runtime::InteropServices::Out]unsigned short% maxvoltage)
- void [SetUseBubble](#) (unsigned short index, unsigned int usebubble)
- unsigned int [GetUseBubble](#) (unsigned short index)
- unsigned short [GetSupplyVoltage](#) ()
- unsigned short [GetAnalogVoltage](#) (unsigned short index)
- unsigned short [GetDigitalIn](#) (unsigned short index)
- unsigned short [GetBubbleState](#) ()

Additional Inherited Members**10.85.1 Constructor & Destructor Documentation**

10.85.1.1 CPPS_FunctionNet() [1/2] [CPPS_FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] cPPS_FunctionPointerContainer)

10.85.1.2 CPPS_FunctionNet() [2/2] [CPPS_FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.85.2 Member Function Documentation

10.85.2.1 GetAnalogVoltage() unsigned short GetAnalogVoltage (
 unsigned short *index*)

10.85.2.2 GetAnalogVoltages() void GetAnalogVoltages (
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned short% *minvoltage*,
 [System::Runtime::InteropServices::Out] unsigned short% *maxvoltage*)

10.85.2.3 GetBubbleState() unsigned short GetBubbleState ()

10.85.2.4 GetDigitalIn() unsigned short GetDigitalIn (
 unsigned short *index*)

10.85.2.5 GetPumpCouple() unsigned int GetPumpCouple ()

10.85.2.6 GetPumpEnableSpeedRatio() unsigned int GetPumpEnableSpeedRatio ()

10.85.2.7 GetPumpFastOnOff() unsigned int GetPumpFastOnOff (
 unsigned short *index*)

10.85.2.8 GetPumpFastSpeed() short GetPumpFastSpeed (
 unsigned short *index*)

10.85.2.9 GetPumpFunctionSpeeds() void GetPumpFunctionSpeeds (
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] short% *offspeed*,
 [System::Runtime::InteropServices::Out] short% *onspeed*)

10.85.2.10 GetPumpManualOnOff() unsigned int GetPumpManualOnOff (
 unsigned short *index*)

10.85.2.11 GetPumpMaxSpeed() unsigned short GetPumpMaxSpeed (
 unsigned short *index*)

10.85.2.12 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType (
 unsigned short *index*)

10.85.2.13 GetPumpSpeedRatio() int GetPumpSpeedRatio ()

10.85.2.14 GetPumpSpeedUnit() int GetPumpSpeedUnit (
 unsigned short *index*)

10.85.2.15 GetSupplyVoltage() unsigned short GetSupplyVoltage ()

10.85.2.16 GetUseBubble() unsigned int GetUseBubble (
 unsigned short *index*)

10.85.2.17 SetAnalogVoltages() void SetAnalogVoltages (
 unsigned short *index*,
 unsigned short *minvoltage*,
 unsigned short *maxvoltage*)

10.85.2.18 SetPumpCouple() void SetPumpCouple (
 unsigned int *i*)

10.85.2.19 SetPumpEnableSpeedRatio() void SetPumpEnableSpeedRatio (
 unsigned int *enable*)

10.85.2.20 SetPumpFastOnOff() void SetPumpFastOnOff (
 unsigned short *index*,
 unsigned int *onoff*)

10.85.2.21 SetPumpFastSpeed() void SetPumpFastSpeed (
 unsigned short *index*,
 short *fastspeed*)

10.85.2.22 SetPumpFunctionSpeeds() void SetPumpFunctionSpeeds (
 unsigned short *index*,
 short *offspeed*,
 short *onspeed*)

10.85.2.23 SetPumpManualOnOff() void SetPumpManualOnOff (
 unsigned short *index*,
 unsigned int *onoff*)

10.85.2.24 SetPumpMaxSpeed() void SetPumpMaxSpeed (
 unsigned short *index*,
 unsigned short *maxspeed*)

10.85.2.25 SetPumpModeType() void SetPumpModeType (
 unsigned short *index*,
 PP_Pump_Mode_Type_EnumNet *type*)

10.85.2.26 SetPumpSpeedRatio() void SetPumpSpeedRatio (
 int *ratio*)


10.85.2.27 SetPumpSpeedUnit() void SetPumpSpeedUnit (
 unsigned short *index*,
 int *speedunit*)

10.85.2.28 SetUseBubble() `void SetUseBubble (`
 `unsigned short index,`
 `unsigned int usebubble)`

10.86 CPPSDeviceNet Class Reference

CPPS4plus1DeviceNet is the to control the MCS HLA device

Inheritance diagram for CPPSDeviceNet:

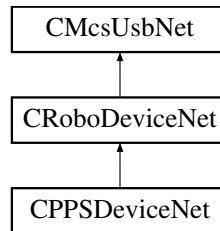

Public Member Functions

- [CPPSDeviceNet](#) (void)

Additional Inherited Members

10.86.1 Detailed Description

CPPS4plus1DeviceNet is the to control the MCS HLA device


10.86.2 Constructor & Destructor Documentation

10.86.2.1 CPPSDeviceNet() `CPPSDeviceNet (`
 `void)`

10.87 CProgramPressureCurveNet Class Reference

[CProgramPressureCurveNet](#) is the class to program pressure curves

Inheritance diagram for CProgramPressureCurveNet:

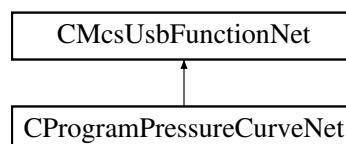

Public Member Functions

- [CProgramPressureCurveNet](#) ([CMcsUsbNet](#)^ mcsusb)
Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.
- virtual [~CProgramPressureCurveNet](#) (void)
- [ICProgramPressureCurveNet](#) (void)
- void [Program](#) (unsigned char busnumber, unsigned char busaddress, int32_t channel, array< int32_t >^ pressures, array< int32_t >^ steps, array< int16_t >^ durations)
- void [SetRepeats](#) (unsigned char busnumber, unsigned char busaddress, int32_t channel, uint32_t repeats)
- unsigned int [GetRepeats](#) (unsigned char busnumber, unsigned char busaddress, int32_t channel)

Additional Inherited Members

10.87.1 Detailed Description

[CProgramPressureCurveNet](#) is the class to program pressure curves

10.87.2 Constructor & Destructor Documentation

10.87.2.1 CProgramPressureCurveNet() [CProgramPressureCurveNet](#) ([CMcsUsbNet](#)^ mcsusb)

Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.

10.87.2.2 ~CProgramPressureCurveNet() virtual [~CProgramPressureCurveNet](#) (void) [virtual]

10.87.2.3 "!ICProgramPressureCurveNet() [!CProgramPressureCurveNet](#) (void)

10.87.3 Member Function Documentation

10.87.3.1 GetRepeats() unsigned int [GetRepeats](#) (unsigned char busnumber, unsigned char busaddress, int32_t channel)


10.87.3.2 Program() void Program (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int32_t *channel*,
 array< int32_t >^ *pressures*,
 array< int32_t >^ *steps*,
 array< int16_t >^ *durations*)

10.87.3.3 SetRepeats() void SetRepeats (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int32_t *channel*,
 uint32_t *repeats*)

10.88 CPulseGeneratorFunctionNet Class Reference

[CPulseGeneratorFunctionNet](#) is the class to control the pulse generator for video tracking

Inheritance diagram for CPulseGeneratorFunctionNet:

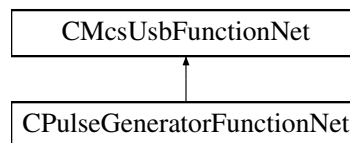

Public Member Functions

- [CPulseGeneratorFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pPulseGeneratorFunctionPointerContainer)
Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.
- [CPulseGeneratorFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual ~[CPulseGeneratorFunctionNet](#) ()
- ![CPulseGeneratorFunctionNet](#) ()
- int32_t [GetPeriod](#) (int32_t generator_number)
Reads the generator period
- void [SetPeriod](#) (int32_t generator_number, int32_t period_in_samples)
Writes the generator period
- int32_t [GetPulseLength](#) (int32_t generator_number)
Reads the generator pulse length
- void [SetPulseLength](#) (int32_t generator_number, int32_t pulselength_in_10us)
Writes the generator pulse length
- void [GetModeSelect](#) (int32_t generator_number, [System::Runtime::InteropServices::Out]PulseGeneratorModeEnumNet% mode, [System::Runtime::InteropServices::Out]int32_t% digitalchannel)
Reads the generator mode
- void [SetModeSelect](#) (int32_t generator_number, PulseGeneratorModeEnumNet mode, int32_t digitalchannel)
Writes the generator mode

Additional Inherited Members

10.88.1 Detailed Description

[CPulseGeneratorFunctionNet](#) is the class to control the pulse generator for video tracking

10.88.2 Constructor & Destructor Documentation

10.88.2.1 CPulseGeneratorFunctionNet() [1/2] [CPulseGeneratorFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *pPulseGeneratorFunctionPointerContainer*)

Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.

10.88.2.2 CPulseGeneratorFunctionNet() [2/2] [CPulseGeneratorFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

10.88.2.3 ~CPulseGeneratorFunctionNet() [virtual ~CPulseGeneratorFunctionNet](#) () [virtual]

10.88.2.4 "!CPulseGeneratorFunctionNet() [!CPulseGeneratorFunctionNet](#) ()

10.88.3 Member Function Documentation

10.88.3.1 GetModeSelect() [void GetModeSelect](#) (
[int32_t](#) *generator_number*,
[\[System::Runtime::InteropServices::Out\]](#) [PulseGenerator_Mode_EnumNet](#)% *mode*,
[\[System::Runtime::InteropServices::Out\]](#) [int32_t](#)% *digitalchannel*)

Reads the generator mode

Parameters

<i>generator_number</i>	The generator number
<i>mode</i>	The generator mode
<i>digitalchannel</i>	The digital in channel used as gate

10.88.3.2 GetPeriod() `int32_t GetPeriod (`
`int32_t generator_number)`

Reads the generator period

Parameters

<i>generator_number</i>	The generator number
-------------------------	----------------------

Returns

The period

10.88.3.3 GetPulseLength() `int32_t GetPulseLength (`
`int32_t generator_number)`

Reads the generator pulse length

Parameters

<i>generator_number</i>	The generator number
-------------------------	----------------------

Returns

The pulse length

10.88.3.4 SetModeSelect() `void SetModeSelect (`
`int32_t generator_number,`
`PulseGenerator_Mode_EnumNet mode,`
`int32_t digitalchannel)`

Writes the generator mode

Parameters

<i>generator_number</i>	The generator number
<i>mode</i>	The generator mode
<i>digitalchannel</i>	The digital in channel used as gate

10.88.3.5 SetPeriod() `void SetPeriod (`

```
int32_t generator_number,
int32_t period_in_samples )
```

Writes the generator period

Parameters

<i>generator_number</i>	The generator number
<i>period_in_samples</i>	The period

10.88.3.6 SetPulseLength() void SetPulseLength (
int32_t generator_number,
int32_t pulselength_in_10us)


Writes the generator pulse length

Parameters

<i>generator_number</i>	The generator number
<i>pulselength_in_10us</i>	The pulse length

10.89 CRadioControlledDevicesNet Class Reference

Inheritance diagram for CRadioControlledDevicesNet:

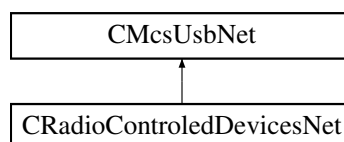

Public Member Functions

- [CRadioControlledDevicesNet](#) (void)
- bool [HasRadioControl](#) ()
- array< unsigned short > ^ [GetDeviceNames](#) ()
- void [ConnectDevice](#) (unsigned short sn)
- void [DisConnectDevice](#) ()
- bool [StillConnected](#) ()
- void [SetFrequency](#) (unsigned short frequency)
- unsigned short [GetFrequency](#) ()

Protected Member Functions

- [CRadioControlledDevicesNet](#) (CRadioControlledDevices *pRadioControlled)

Additional Inherited Members

10.89.1 Constructor & Destructor Documentation

10.89.1.1 CRadioControlledDevicesNet() [1/2] `CRadioControlledDevicesNet (void)`

10.89.1.2 CRadioControlledDevicesNet() [2/2] `CRadioControlledDevicesNet (CRadioControlledDevices * pRadioControlled) [protected]`

10.89.2 Member Function Documentation

10.89.2.1 ConnectDevice() `void ConnectDevice (unsigned short sn)`

10.89.2.2 DisconnectDevice() `void DisconnectDevice ()`

10.89.2.3 GetDeviceNames() `array<unsigned short> ^ GetDeviceNames ()`

10.89.2.4 GetFrequency() `unsigned short GetFrequency ()`

10.89.2.5 HasRadioControl() `bool HasRadioControl ()`

10.89.2.6 SetFrequency() `void SetFrequency (unsigned short frequency)`

10.89.2.7 StillConnected() `bool StillConnected ()`

10.90 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference

Public Member Functions

- [CRegionOfInterestRect](#) (int left, int top, int right, int bottom)
- [CRegionOfInterestRect](#) ^ [DeepCopy](#) ()

Public Attributes

- int [m_Left](#)
- int [m_Top](#)
- int [m_Right](#)
- int [m_Bottom](#)

10.90.1 Constructor & Destructor Documentation

10.90.1.1 [CRegionOfInterestRect](#)() [CRegionOfInterestRect](#) (
 int *left*,
 int *top*,
 int *right*,
 int *bottom*)

10.90.2 Member Function Documentation

10.90.2.1 [DeepCopy](#)() [CRegionOfInterestRect](#) ^ [DeepCopy](#) ()

10.90.3 Member Data Documentation

10.90.3.1 [m_Bottom](#) int m_Bottom


10.90.3.2 [m_Left](#) int m_Left

10.90.3.3 [m_Right](#) int m_Right

10.90.3.4 m_Top `int m_Top`

10.91 CRetinaLedDeviceNet Class Reference

Inheritance diagram for CRetinaLedDeviceNet:

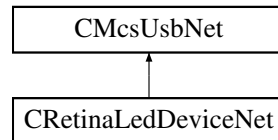

Public Member Functions

- [CRetinaLedDeviceNet](#) ()
- [~CRetinaLedDeviceNet](#) ()
- unsigned int [SetTrigger](#) (int enable)
- unsigned int [SetLED](#) (unsigned long long pattern)
- unsigned int [SetTablepointer](#) (int position)
- unsigned int [GetTablepointer](#) (int % position)
- unsigned int [ClearTable](#) ()
- unsigned int [AddTableEntry](#) (unsigned long long pattern)
- unsigned int [AddLoopEntry](#) (unsigned short repeats, unsigned short steps_back)
- unsigned int [SetRepeat](#) (int repeat)
- unsigned int [SetLumi](#) (int lumi)
- unsigned int [SetPersistency](#) (unsigned int persistency)

Additional Inherited Members

10.91.1 Constructor & Destructor Documentation

10.91.1.1 CRetinaLedDeviceNet() [CRetinaLedDeviceNet](#) ()

10.91.1.2 ~CRetinaLedDeviceNet() [~CRetinaLedDeviceNet](#) ()

10.91.2 Member Function Documentation

10.91.2.1 AddLoopEntry() unsigned int AddLoopEntry (
 unsigned short *repeats*,
 unsigned short *steps_back*)

10.91.2.2 AddTableEntry() unsigned int AddTableEntry (
 unsigned long long *pattern*)

10.91.2.3 ClearTable() unsigned int ClearTable ()

10.91.2.4 GetTablepointer() unsigned int GetTablepointer (
 int % *position*)

10.91.2.5 SetLED() unsigned int SetLED (
 unsigned long long *pattern*)

10.91.2.6 SetLumi() unsigned int SetLumi (
 int *lumi*)

10.91.2.7 SetPersistency() unsigned int SetPersistency (
 unsigned int *persistency*)

10.91.2.8 SetRepeat() unsigned int SetRepeat (
 int *repeat*)


10.91.2.9 SetTablepointer() unsigned int SetTablepointer (
 int *position*)

10.91.2.10 SetTrigger() unsigned int SetTrigger (
 int *enable*)

10.92 CRFFunctionNet Class Reference

[CRFFunctionNet](#) is the class to control RF devices

Inheritance diagram for CRFFunctionNet:

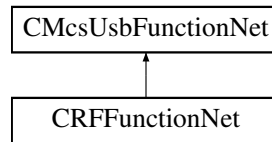

Public Member Functions

- [CRFFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pRFFunctionPointerContainer)
 - Initializes a new instance of the [CRFFunctionNet](#) class.*
- [CRFFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CRFFunctionNet](#) ()
- [!CRFFunctionNet](#) ()
- uint32_t [GetBaseFrequency](#) (CFirmwareDestinationNet destination)
 - gets the base advertise frequency*
- void [SetBaseFrequency](#) (CFirmwareDestinationNet destination, uint32_t frequency)
 - sets the base advertise frequency*
- uint32_t [GetWorkingFrequency](#) ()
 - gets the working frequency*
- void [SetWorkingFrequency](#) (uint32_t frequency)
 - sets the working frequency*
- array< uint32_t > [^] [GetAvailableDeviceListEx](#) (int list_Length)
 - get a list of available devices*
- array< uint32_t > [^] [GetAvailableDeviceList](#) ()
 - get a list of available devices*
- array< uint32_t > [^] [GetAvailableStateListEx](#) (int list_Length)
 - get a list of the states of the available devices*
- array< uint32_t > [^] [GetAvailableStateList](#) ()
 - get a list of the states of the available devices*
- void [Connect](#) (uint32_t sn)
 - connect to a RF device, use 0 to disconnect*
- uint32_t [GetConnectedDevice](#) ()
 - get connect RF device, 0 = no device connected*
- uint32_t [GetState](#) ()
 - get connection state*
- void [SetTestMode](#) (uint32_t mode)
 - set test mode*
- uint32_t [GetTestMode](#) ()
 - gets test mode*

Additional Inherited Members

10.92.1 Detailed Description

[CRFFunctionNet](#) is the class to control RF devices

10.92.2 Constructor & Destructor Documentation

10.92.2.1 CRFFunctionNet() [1/2] `CRFFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pRFFunctionPointerContainer)`

Initializes a new instance of the `CRFFunctionNet` class.

10.92.2.2 CRFFunctionNet() [2/2] `CRFFunctionNet (CMcsUsbNet^ mcsusb)`

10.92.2.3 ~CRFFunctionNet() `virtual ~CRFFunctionNet () [virtual]`

10.92.2.4 !CRFFunctionNet() `!CRFFunctionNet ()`

10.92.3 Member Function Documentation

10.92.3.1 Connect() `void Connect (uint32_t sn)`

connect to a RF device, use 0 to disconnect

Parameters

<i>sn</i>	
-----------	--

10.92.3.2 GetAvailableDeviceList() `array<uint32_t> ^ GetAvailableDeviceList ()`

get a list of available devices

Returns

10.92.3.3 GetAvailableDeviceListEx() `array<uint32_t> ^ GetAvailableDeviceListEx (`
`int list_Length)`

get a list of available devices

Parameters

<i>list_Length</i>	The maximal length of list.
--------------------	-----------------------------

Returns

10.92.3.4 GetAvailableStateList() `array<uint32_t> ^ GetAvailableStateList ()`

get a list of the states of the available devices

Returns

10.92.3.5 GetAvailableStateListEx() `array<uint32_t> ^ GetAvailableStateListEx (`
`int list_Length)`

get a list of the states of the available devices

Parameters

<i>list_Length</i>	The maximal length of list.
--------------------	-----------------------------

Returns

10.92.3.6 GetBaseFrequency() `uint32_t GetBaseFrequency (`
`CFirmwareDestinationNet destination)`

gets the base advertise frequency

Parameters

<i>destination</i>	
--------------------	--

Returns

10.92.3.7 GetConnectedDevice() `uint32_t GetConnectedDevice ()`

get connect RF device, 0 = no device connected

Returns

10.92.3.8 GetState() `uint32_t GetState ()`

get connection state

Returns

10.92.3.9 GetTestMode() `uint32_t GetTestMode ()`

gets test mode

Returns

10.92.3.10 GetWorkingFrequency() `uint32_t GetWorkingFrequency ()`

gets the working frequency

Returns

10.92.3.11 SetBaseFrequency() `void SetBaseFrequency (`
 `CFirmwareDestinationNet destination,`
 `uint32_t frequency)`

gets the base advertise frequency

Parameters

<i>destination</i>	
<i>frequency</i>	

10.92.3.12 SetTestMode() void SetTestMode (
 uint32_t mode)

set test mode

Parameters

<i>mode</i>	
-------------	--

10.92.3.13 SetWorkingFrequency() void SetWorkingFrequency (
 uint32_t frequency)


sets the working frequency

Parameters

<i>frequency</i>	
------------------	--

10.93 CRobo_FYIProgram_FunctionNet Class Reference

Inheritance diagram for CRobo_FYIProgram_FunctionNet:

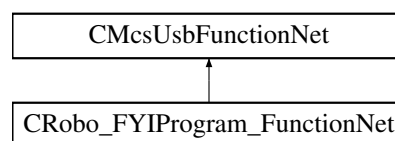

Public Member Functions

- [CRobo_FYIProgram_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] robo_FY↔ IProgram_FunctionPointerContainer)
- [CRobo_FYIProgram_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetValve1](#) (unsigned char index, unsigned int valve1)
- unsigned int [GetValve1](#) (unsigned char index)
- void [SetValve2](#) (unsigned char index, unsigned int valve2)
- unsigned int [GetValve2](#) (unsigned char index)
- void [SetLength](#) (unsigned char index, int length)
- int [GetLength](#) (unsigned char index)
- void [Start](#) ()
- int [GetState](#) ()

Additional Inherited Members

10.93.1 Constructor & Destructor Documentation

10.93.1.1 CRobo_FYIProgram_FunctionNet() [1/2] CRobo_FYIProgram_FunctionNet (
 CMcsUsbNet^ *mcsusb*,
 CMcsUsbFunctionPointerContainer^ *robo_FYIProgram_FunctionPointerContainer*)

10.93.1.2 CRobo_FYIProgram_FunctionNet() [2/2] CRobo_FYIProgram_FunctionNet (
 CMcsUsbNet^ *mcsusb*)

10.93.2 Member Function Documentation

10.93.2.1 GetLength() int GetLength (
 unsigned char *index*)

10.93.2.2 GetState() int GetState ()

10.93.2.3 GetValve1() unsigned int GetValve1 (
 unsigned char *index*)

10.93.2.4 GetValve2() unsigned int GetValve2 (
 unsigned char *index*)

10.93.2.5 SetLength() void SetLength (
 unsigned char *index*,
 int *length*)


10.93.2.6 SetValve1() void SetValve1 (
 unsigned char *index*,
 unsigned int *valve1*)

10.93.2.7 SetValve2() void SetValve2 (
 unsigned char *index*,
 unsigned int *valve2*)

10.93.2.8 Start() void Start ()

10.94 CRobo_FYITemp_FunctionNet Class Reference

Inheritance diagram for CRobo_FYITemp_FunctionNet:

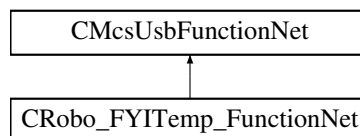

Public Member Functions

- [CRobo_FYITemp_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetRegulatorOnOff](#) (unsigned char *index*, int *onoff*)
- int [GetRegulatorOnOff](#) (unsigned char *index*)
- void [SetSollTemp](#) (unsigned char *index*, int *temp*)
- int [GetSollTemp](#) (unsigned char *index*)
- void [SetPCoeff](#) (unsigned char *index*, int *pcoeff*)
- int [GetPCoeff](#) (unsigned char *index*)
- void [SetICoeff](#) (unsigned char *index*, int *icoeff*)
- int [GetICoeff](#) (unsigned char *index*)
- void [SetMaxPower](#) (unsigned char *index*, int *power*)
- int [GetMaxPower](#) (unsigned char *index*)

Additional Inherited Members

10.94.1 Constructor & Destructor Documentation

10.94.1.1 CRobo_FYITemp_FunctionNet() [CRobo_FYITemp_FunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

10.94.2 Member Function Documentation

10.94.2.1 GetICoeff() `int GetICoeff (unsigned char index)`

10.94.2.2 GetMaxPower() `int GetMaxPower (unsigned char index)`

10.94.2.3 GetPCoeff() `int GetPCoeff (unsigned char index)`

10.94.2.4 GetRegulatorOnOff() `int GetRegulatorOnOff (unsigned char index)`

10.94.2.5 GetSoilTemp() `int GetSoilTemp (unsigned char index)`

10.94.2.6 SetICoeff() `void SetICoeff (unsigned char index,
int icoeff)`

10.94.2.7 SetMaxPower() `void SetMaxPower (unsigned char index,
int power)`


10.94.2.8 SetPCoeff() `void SetPCoeff (unsigned char index,
int pcoeff)`

10.94.2.9 SetRegulatorOnOff() void SetRegulatorOnOff (
 unsigned char *index*,
 int *onoff*)

10.94.2.10 SetSollTemp() void SetSollTemp (
 unsigned char *index*,
 int *temp*)

10.95 CRoboDacqNet Class Reference

Inheritance diagram for CRoboDacqNet:

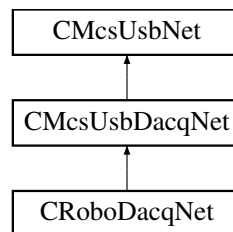

Public Member Functions

- [CRoboDacqNet](#) (void)
- [CRoboDacqNet](#) ([CRoboDeviceNet](#)[^] robodevice)
- void [RunTable](#) ()
- void [RunTable](#) (int timeout)
- void [StopTable](#) ()
- void [StopTable](#) (int timeout)
- void [CancelTableLoop](#) ()
- void [CancelTableLoopAndStopTable](#) ()
- void [SetConfigurationBit](#) (unsigned short bit, bool value)
- void [SetConfigurationBitSupply](#) (bool value)
- void [SetConfigurationBitRelais](#) (bool value)
- void [SetConfigurationBitStream](#) (bool value)
- void [SetConfigurationBitAxc](#) (bool value)
- void [SetConfigurationBitCC_Gen](#) (bool value)
- void [SetConfigurationBitCV_Gen](#) (bool value)
- void [SetConfigurationBitRC_Gen](#) (bool value)
- void [SetConfigurationBitRV_Gen](#) (bool value)
- void [SetConfigurationBitBlu_Led](#) (bool value)
- void [SetConfigurationBitRed_Led](#) (bool value)
- void [SetConfigurationBitBlu_LedToggleSlow](#) (bool value)
- void [SetConfigurationBitRed_LedToggleSlow](#) (bool value)
- void [SetConfigurationBitBlu_LedToggleFast](#) (bool value)
- void [SetConfigurationBitRed_LedToggleFast](#) (bool value)
- void [SetConfigurationBitRed_LedSaturation](#) (bool value)
- void [SetSimulation](#) (unsigned int enable)
- void [SetUClamp](#) (int uClamp)
- void [SetIClamp](#) (int iClamp)

- void [SetPGain](#) (int pGain)
- void [SetIGain](#) (int iGain)
- void [SetFilter](#) (int filter)
- void [SetUJOffset](#) (int UJOffset)
- void [SetUCOffset](#) (int UCOffset)
- void [SetICOffset](#) (int ICOffset)
- void [SetCrossTalkOffset](#) (int CrossTalk)
- void [SetXGain](#) (int xGain)
- void [SetCrossTalkOptimum](#) (int cxOptimum)
- void [SetRecordingNumber](#) (unsigned int recordingNumber)
- void [ClampAmpRestart](#) ()
- void [DoRamp](#) (int startValue, int endValue, int duration, int mode)
- unsigned int [GetClampAmpSerialNumber](#) ()
- unsigned int [GetConfigurationBits](#) ()
- bool [GetConfigurationBit](#) (unsigned short bit)
- bool [GetConfigurationBitSupply](#) ()
- bool [GetConfigurationBitRelais](#) ()
- bool [GetConfigurationBitStream](#) ()
- bool [GetConfigurationBitAxc](#) ()
- bool [GetConfigurationBitCC_Gen](#) ()
- bool [GetConfigurationBitCV_Gen](#) ()
- bool [GetConfigurationBitRC_Gen](#) ()
- bool [GetConfigurationBitRV_Gen](#) ()
- bool [GetConfigurationBitBlu_Led](#) ()
- bool [GetConfigurationBitRed_Led](#) ()
- bool [GetConfigurationBitBlu_LedToggleSlow](#) ()
- bool [GetConfigurationBitRed_LedToggleSlow](#) ()
- bool [GetConfigurationBitBlu_LedToggleFast](#) ()
- bool [GetConfigurationBitRed_LedToggleFast](#) ()
- bool [GetConfigurationBitRed_LedSaturation](#) ()
- unsigned int [GetSimulation](#) ()
- int [GetUClamp](#) ()
- int [GetIClamp](#) ()
- int [GetPGain](#) ()
- int [GetIGain](#) ()
- int [GetFilter](#) ()
- int [GetUJOffset](#) ()
- int [GetUCOffset](#) ()
- int [GetICOffset](#) ()
- int [GetCrossTalkOffset](#) ()
- int [GetXGain](#) ()
- int [GetCrossTalkOptimum](#) ()
- unsigned int [GetRecordingNumber](#) ()
- int [GetResistanceC](#) ()
- int [GetResistanceV](#) ()
- int [GetCapacityC](#) ()
- int [GetCapacityV](#) ()
- int [GetCapacityX](#) ()
- int [GetUV](#) ()
- int [GetUC](#) ()
- int [GetIC](#) ()
- int [GetNUV_MS](#) ()
- int [GetNUC_MS](#) ()
- int [GetNIC_MS](#) ()
- void [SetAllDigout](#) (uint32_t value)

- uint32_t [GetAllDigout](#) ()
- void [SetCommand](#) (unsigned char command, int value)
- int [GetCommand](#) (unsigned char command)
- void [SetDigout](#) (uint16_t index, bool enable)
- bool [GetDigout](#) (uint16_t index)
- void [TableDefBegin](#) ()
- void [TableDefEnd](#) ()
- void [Table_Wait](#) (unsigned int tableWait)
- void [SetDownsampleFactor](#) (int index, int downsample_factor)
- void [SetFilterCoeffs](#) (int index, array< int > ^ coeffs)
- int [GetDownsampleFactor](#) (int index)
- array< int > ^ [GetFilterCoeffs](#) (int index)
- void [Emu_SetElectrodeResists](#) (int emuElectrodeResist)
- void [Emu_SetCellResists](#) (int emuCellResist)
- void [Emu_SetCellCapacity](#) (int emuCellCapacity)
- void [Emu_SetCellPotential](#) (int emuCellPotential)
- void [Emu_SetNoise](#) (int emuNoise)
- int [Emu_GetElectrodeResists](#) ()
- int [Emu_GetCellResists](#) ()
- int [Emu_GetCellCapacity](#) ()
- int [Emu_GetCellPotential](#) ()
- int [Emu_GetNoise](#) ()
- void [SetDisplayText](#) (int index, String ^ displayText)
- void [SetScreen](#) (int screen)
- void [UpdateDisplay](#) ()
- String ^ [GetDisplayText](#) (int index)
- int [GetScreen](#) ()
- int [GetUpdateDisplay](#) ()

Additional Inherited Members

10.95.1 Constructor & Destructor Documentation

10.95.1.1 CRoboDacqNet() [1/2] [CRoboDacqNet](#) (
 void)

10.95.1.2 CRoboDacqNet() [2/2] [CRoboDacqNet](#) (
 [CRoboDeviceNet](#) ^ robodevice)

10.95.2 Member Function Documentation

10.95.2.1 CancelTableLoop() void CancelTableLoop ()

10.95.2.2 CancelTableLoopAndStopTable() `void CancelTableLoopAndStopTable ()`

10.95.2.3 ClampAmpRestart() `void ClampAmpRestart ()`

10.95.2.4 DoRamp() `void DoRamp (`
 `int startValue,`
 `int endValue,`
 `int duration,`
 `int mode)`

10.95.2.5 Emu_GetCellCapacity() `int Emu_GetCellCapacity ()`

10.95.2.6 Emu_GetCellPotential() `int Emu_GetCellPotential ()`

10.95.2.7 Emu_GetCellResists() `int Emu_GetCellResists ()`

10.95.2.8 Emu_GetElectrodeResists() `int Emu_GetElectrodeResists ()`

10.95.2.9 Emu_GetNoise() `int Emu_GetNoise ()`

10.95.2.10 Emu_SetCellCapacity() `void Emu_SetCellCapacity (`
 `int emuCellCapacity)`

10.95.2.11 Emu_SetCellPotential() `void Emu_SetCellPotential (`
 `int emuCellPotential)`

10.95.2.12 Emu_SetCellResists() void Emu_SetCellResists (
int *emuCellResist*)

10.95.2.13 Emu_SetElectrodeResists() void Emu_SetElectrodeResists (
int *emuElectrodeResist*)

10.95.2.14 Emu_SetNoise() void Emu_SetNoise (
int *emuNoise*)

10.95.2.15 GetAllDigout() uint32_t GetAllDigout ()

10.95.2.16 GetCapacityC() int GetCapacityC ()

10.95.2.17 GetCapacityV() int GetCapacityV ()

10.95.2.18 GetCapacityX() int GetCapacityX ()

10.95.2.19 GetClampAmpSerialNumber() unsigned int GetClampAmpSerialNumber ()

10.95.2.20 GetCommand() int GetCommand (
unsigned char *command*)

10.95.2.21 GetConfigurationBit() bool GetConfigurationBit (
unsigned short *bit*)

10.95.2.22 GetConfigurationBitAxc() `bool GetConfigurationBitAxc ()`

10.95.2.23 GetConfigurationBitBlu_Led() `bool GetConfigurationBitBlu_Led ()`

10.95.2.24 GetConfigurationBitBlu_LedToggleFast() `bool GetConfigurationBitBlu_LedToggleFast ()`

10.95.2.25 GetConfigurationBitBlu_LedToggleSlow() `bool GetConfigurationBitBlu_LedToggleSlow ()`

10.95.2.26 GetConfigurationBitCC_Gen() `bool GetConfigurationBitCC_Gen ()`

10.95.2.27 GetConfigurationBitCV_Gen() `bool GetConfigurationBitCV_Gen ()`

10.95.2.28 GetConfigurationBitRC_Gen() `bool GetConfigurationBitRC_Gen ()`

10.95.2.29 GetConfigurationBitRed_Led() `bool GetConfigurationBitRed_Led ()`

10.95.2.30 GetConfigurationBitRed_LedSaturation() `bool GetConfigurationBitRed_LedSaturation ()`

10.95.2.31 GetConfigurationBitRed_LedToggleFast() `bool GetConfigurationBitRed_LedToggleFast ()`

10.95.2.32 GetConfigurationBitRed_LedToggleSlow() `bool GetConfigurationBitRed_LedToggleSlow ()`

10.95.2.33 GetConfigurationBitRelais() `bool GetConfigurationBitRelais ()`

10.95.2.34 GetConfigurationBitRV_Gen() `bool GetConfigurationBitRV_Gen ()`

10.95.2.35 GetConfigurationBits() `unsigned int GetConfigurationBits ()`

10.95.2.36 GetConfigurationBitStream() `bool GetConfigurationBitStream ()`

10.95.2.37 GetConfigurationBitSupply() `bool GetConfigurationBitSupply ()`

10.95.2.38 GetCrossTalkOffset() `int GetCrossTalkOffset ()`

10.95.2.39 GetCrossTalkOptimum() `int GetCrossTalkOptimum ()`

10.95.2.40 GetDigout() `bool GetDigout (`
`uint16_t index)`

10.95.2.41 GetDisplayText() `String ^ GetDisplayText (`
`int index)`

10.95.2.42 GetDownsampleFactor() `int GetDownsampleFactor (`
`int index)`

10.95.2.43 GetFilter() `int GetFilter ()`

10.95.2.44 GetFilterCoeffs() `array<int> ^ GetFilterCoeffs (`
`int index)`

10.95.2.45 GetIC() `int GetIC ()`

10.95.2.46 GetIClamp() `int GetIClamp ()`

10.95.2.47 GetICOffset() `int GetICOffset ()`

10.95.2.48 GetIGain() `int GetIGain ()`

10.95.2.49 GetNIC_MS() `int GetNIC_MS ()`

10.95.2.50 GetNUC_MS() `int GetNUC_MS ()`

10.95.2.51 GetNUV_MS() `int GetNUV_MS ()`

10.95.2.52 GetPGain() `int GetPGain ()`

10.95.2.53 GetRecordingNumber() `unsigned int GetRecordingNumber ()`

10.95.2.54 GetResistanceC() `int GetResistanceC ()`

10.95.2.55 GetResistanceV() `int GetResistanceV ()`

10.95.2.56 GetScreen() `int GetScreen ()`

10.95.2.57 GetSimulation() `unsigned int GetSimulation ()`

10.95.2.58 GetUC() `int GetUC ()`

10.95.2.59 GetUClamp() `int GetUClamp ()`

10.95.2.60 GetUCOffset() `int GetUCOffset ()`

10.95.2.61 GetUpdateDisplay() `int GetUpdateDisplay ()`

10.95.2.62 GetUV() `int GetUV ()`

10.95.2.63 GetUVOffset() `int GetUVOffset ()`

10.95.2.64 GetXGain() `int GetXGain ()`

10.95.2.65 RunTable() `[1/2] void RunTable ()`

10.95.2.66 RunTable() [2/2] void RunTable (
int *timeout*)

10.95.2.67 SetAllDigout() void SetAllDigout (
uint32_t *value*)

10.95.2.68 SetCommand() void SetCommand (
unsigned char *command*,
int *value*)

10.95.2.69 SetConfigurationBit() void SetConfigurationBit (
unsigned short *bit*,
bool *value*)

10.95.2.70 SetConfigurationBitAxc() void SetConfigurationBitAxc (
bool *value*)

10.95.2.71 SetConfigurationBitBlu_Led() void SetConfigurationBitBlu_Led (
bool *value*)

10.95.2.72 SetConfigurationBitBlu_LedToggleFast() void SetConfigurationBitBlu_LedToggleFast (
bool *value*)

10.95.2.73 SetConfigurationBitBlu_LedToggleSlow() void SetConfigurationBitBlu_LedToggleSlow (
bool *value*)

10.95.2.74 SetConfigurationBitCC_Gen() void SetConfigurationBitCC_Gen (
bool *value*)

10.95.2.75 SetConfigurationBitCV_Gen() void SetConfigurationBitCV_Gen (
bool value)

10.95.2.76 SetConfigurationBitRC_Gen() void SetConfigurationBitRC_Gen (
bool value)

10.95.2.77 SetConfigurationBitRed_Led() void SetConfigurationBitRed_Led (
bool value)

10.95.2.78 SetConfigurationBitRed_LedSaturation() void SetConfigurationBitRed_LedSaturation (
bool value)

10.95.2.79 SetConfigurationBitRed_LedToggleFast() void SetConfigurationBitRed_LedToggleFast (
bool value)

10.95.2.80 SetConfigurationBitRed_LedToggleSlow() void SetConfigurationBitRed_LedToggleSlow (
bool value)

10.95.2.81 SetConfigurationBitRelais() void SetConfigurationBitRelais (
bool value)

10.95.2.82 SetConfigurationBitRV_Gen() void SetConfigurationBitRV_Gen (
bool value)

10.95.2.83 SetConfigurationBitStream() void SetConfigurationBitStream (
bool value)

10.95.2.84 SetConfigurationBitSupply() void SetConfigurationBitSupply (
bool value)

10.95.2.85 SetCrossTalkOffset() void SetCrossTalkOffset (
int *CrossTalk*)

10.95.2.86 SetCrossTalkOptimum() void SetCrossTalkOptimum (
int *cxOptimum*)

10.95.2.87 SetDigout() void SetDigout (
uint16_t *index*,
bool *enable*)

10.95.2.88 SetDisplayText() void SetDisplayText (
int *index*,
String^ *displayText*)

10.95.2.89 SetDownsampleFactor() void SetDownsampleFactor (
int *index*,
int *downsample_factor*)

10.95.2.90 SetFilter() void SetFilter (
int *filter*)

10.95.2.91 SetFilterCoeffs() void SetFilterCoeffs (
int *index*,
array< int >^ *coeffs*)

10.95.2.92 SetIClamp() void SetIClamp (
int *iClamp*)

10.95.2.93 SetICOffset() void SetICOffset (
int *ICOffset*)

10.95.2.94 SetIGain() void SetIGain (
int *iGain*)

10.95.2.95 SetPGain() void SetPGain (
int *pGain*)

10.95.2.96 SetRecordingNumber() void SetRecordingNumber (
unsigned int *recordingNumber*)

10.95.2.97 SetScreen() void SetScreen (
int *screen*)

10.95.2.98 SetSimulation() void SetSimulation (
unsigned int *enable*)

10.95.2.99 SetUClamp() void SetUClamp (
int *uClamp*)

10.95.2.100 SetUCOffset() void SetUCOffset (
int *UCOffset*)

10.95.2.101 SetUVOffset() void SetUVOffset (
int *UVOffset*)

10.95.2.102 SetXGain() void SetXGain (
int *xGain*)

10.95.2.103 StopTable() [1/2] void StopTable ()

10.95.2.104 StopTable() [2/2] `void StopTable (`
`int timeout)`

10.95.2.105 Table_Wait() `void Table_Wait (`
`unsigned int tableWait)`

10.95.2.106 TableDefBegin() `void TableDefBegin ()`


10.95.2.107 TableDefEnd() `void TableDefEnd ()`

10.95.2.108 UpdateDisplay() `void UpdateDisplay ()`

10.96 CRoboDeviceNet Class Reference

[CRoboDeviceNet](#) is the base class for all Robo platform based devices

Inheritance diagram for CRoboDeviceNet:

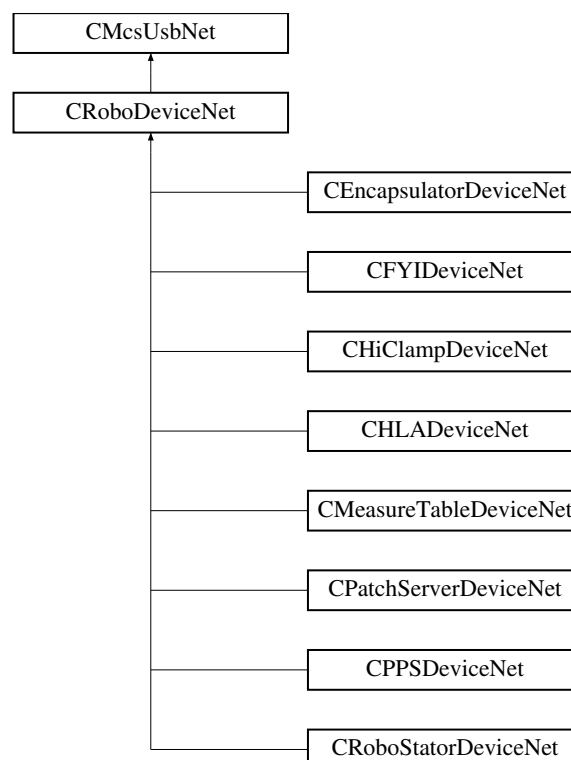

Classes

- class [RoboMainLowLevelCommands](#)

Public Member Functions

- [CRoboDeviceNet](#) (void)
- [~CRoboDeviceNet](#) (void)
- void [SetInMovement](#) ()
Low level command, sets the internal state to "In Movement"
- bool [GetInMovement](#) ()
Low level command, gets the internal state "In Movement"
- uint32_t [GetMovementError](#) ()
Low level command, gets the error of the last movement end
- void [FindReference](#) (unsigned char busaddress, char axes)
- void [FindReference](#) (unsigned char busaddress, char axes, int timeout)
Searches the reference position of the motor
- void [MoveAbs](#) (unsigned char busaddress, char axes, int x, int y)
- void [MoveAbs](#) (unsigned char busaddress, char axes, int x, int y, int timeout)
Moves the motor to the new absolute position
- void [StopMovement](#) (unsigned char busaddress, char axes)
- void [StopMovement](#) (unsigned char busaddress, char axes, int timeout)
Stops the current movement
- void [SetCurrentAndAir](#) (unsigned char busaddress, char axes, unsigned short onoff)
- void [SetCurrentAndAir](#) (unsigned char busaddress, char axes, unsigned short onoff, int timeout)
- void [CancelPoolLoop](#) ()
- void [CancelPoolLoopAndStopMovement](#) ()
- void [GetCurrentPosition](#) (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out]int% x, [System::Runtime::InteropServices::Out]int% y)
Gets the current position of motors
- void [SetAirValve](#) (unsigned int onoff)
- unsigned int [GetAirValve](#) ()
- unsigned int [GetVoltageValves](#) ()
- unsigned int [GetVoltageRs485A](#) ()
- unsigned int [GetVoltageRs485B](#) ()
- unsigned int [GetVoltageAirvalve](#) ()
- unsigned int [GetCurrentAirvalve](#) ()
- unsigned int [GetVoltage12V](#) ()
- unsigned int [GetAirpressure](#) ()
- unsigned int [GetVoltage5V](#) ()
- unsigned int [GetErrorVoltageValves](#) ()
- unsigned int [GetErrorVoltageRs485A](#) ()
- unsigned int [GetErrorVoltageRs485B](#) ()
- unsigned int [GetErrorVoltageAirvalve](#) ()
- unsigned int [GetErrorCurrentAirvalve](#) ()
- unsigned int [GetErrorVoltage12V](#) ()
- unsigned int [GetErrorAirpressure](#) ()
- unsigned int [GetErrorVoltage5V](#) ()
- void [SetVoltageValvesLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageRs485ALimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageRs485BLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageAirvalveLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetCurrentAirvalveLimit](#) (unsigned int lowercurrent, unsigned int uppercurrent)

- void [SetVoltage12VLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetAirpressureLimit](#) (unsigned int lowerpressure, unsigned int upperpressure)
- void [SetVoltage5VLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [GetVoltageValvesLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetVoltageRs485ALimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetVoltageRs485BLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetVoltageAirvalveLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetCurrentAirvalveLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowercurrent, [System::Runtime::InteropServices::Out] unsigned int% uppercurrent)
- void [GetVoltage12VLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetAirpressureLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowerpressure, [System::Runtime::InteropServices::Out] unsigned int% upperpressure)
- void [GetVoltage5VLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [SetMinPressure](#) (int pressure)
- int [GetMinPressure](#) ()

Static Public Attributes

- static const uint32_t [RoboError_Base](#) = (0xA0110000L)
- static const uint32_t [RoboError_UnknownCommand](#) = ((0xA0110000L))
- static const uint32_t [RoboError_Timeout](#) = ((0xA0110000L) | 0x0001)
- static const uint32_t [RoboError_Pressure](#) = ((0xA0110000L) | 0x0002)
- static const uint32_t [RoboError_RangeExceeded](#) = ((0xA0110000L) | 0x0003)
- static const uint32_t [RoboError_CommunicationTimeout](#) = ((0xA0110000L) | 0x0004)
- static const uint32_t [RoboError_AnotherMaster](#) = ((0xA0110000L) | 0x0005)
- static const uint32_t [RoboError_FindReferenceMethod](#) = ((0xA0110000L) | 0x0006)
- static const uint32_t [RoboError_NoSpeedOrAcceleration](#) = ((0xA0110000L) | 0x0007)
- static const uint32_t [RoboError_NoEndSwitch](#) = ((0xA0110000L) | 0x0008)
- static const uint32_t [RoboError_CannotEscapeEndSwitch](#) = ((0xA0110000L) | 0x0009)
- static const uint32_t [RoboError_CommandAlreadyInProgress](#) = ((0xA0110000L) | 0x000A)
- static const uint32_t [RoboError_NoReference](#) = ((0xA0110000L) | 0x000B)
- static const uint32_t [RoboError_OverPressure](#) = ((0xA0110000L) | 0x000C)
- static const uint32_t [RoboError_Phase0OutOfRange](#) = ((0xA0110000L) | 0x000D)
- static const uint32_t [RoboError_PeristalticTimeout](#) = ((0xA0110000L) | 0x000E)
- static const uint32_t [RoboError_GilsonTimeout](#) = ((0xA0110000L) | 0x000F)
- static const uint32_t [RoboError_GilsonWrondID](#) = ((0xA0110000L) | 0x0010)
- static const uint32_t [RoboError_GilsonCommandPending](#) = ((0xA0110000L) | 0x0011)
- static const uint32_t [RoboError_ParameterNotAllowed](#) = ((0xA0110000L) | 0x0012)
- static const uint32_t [RoboError_StateChangeNotPossible](#) = ((0xA0110000L) | 0x0013)
- static const uint32_t [RoboError_CommandNotPossible](#) = ((0xA0110000L) | 0x0014)
- static const uint32_t [RoboError_DacqNotReady](#) = ((0xA0110000L) | 0x0015)
- static const uint32_t [RoboError_NoMoreData](#) = ((0xA0110000L) | 0x0016)
- static const uint32_t [RoboError_McsBus_UnknownCommand](#) = ((0xA0110000L) | 0x003F)
- static const uint32_t [RoboError_DLLMovementTimeout](#) = ((0xA0110000L) | 0x1001)
- static const uint32_t [RoboError_PollLoopCanceled](#) = ((0xA0110000L) | 0x1002)
- static const uint32_t [RoboError_PollLoopCanceledAndStopMovement](#) = ((0xA0110000L) | 0x1003)
- static const byte [McsBus_XY](#) = 1

McsBus address for the xy-plane

- static const byte [McsBus_ZI](#) = 2
McsBus address for the z and i axes
- static const byte [Axis_X](#) = 0
Axis number of x for axis argument
- static const byte [Axis_Y](#) = 1
Axis number of y for axis argument
- static const byte [Axis_Z](#) = 0
Axis number of z for axis argument
- static const byte [Axis_I](#) = 1
Axis number of i for axis argument
- static const char [Axes_X](#) = 1
Bit pattern for x axis for axes argument
- static const char [Axes_Y](#) = 2
Bit pattern for y axis for axes argument
- static const char [Axes_Z](#) = 1
Bit pattern for z axis for axes argument
- static const char [Axes_I](#) = 2
Bit pattern for i axis for axes argument

Properties

- [CMcsBusNet](#)[^] [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]
- [RoboMainLowLevelCommands](#)[^] [RoboMainLowLevelCommand](#) [get]

Events

- [RoboStatusEventDelegate](#)[^] [RoboStatusEvent](#)

Additional Inherited Members

10.96.1 Detailed Description

[CRoboDeviceNet](#) is the base class for all Robo platform based devices

10.96.2 Constructor & Destructor Documentation

10.96.2.1 CRoboDeviceNet() [CRoboDeviceNet](#) (
void)

10.96.2.2 ~CRoboDeviceNet() [~CRoboDeviceNet](#) (
void)

10.96.3 Member Function Documentation

10.96.3.1 CancelPoolLoop() `void CancelPoolLoop ()`

10.96.3.2 CancelPoolLoopAndStopMovement() `void CancelPoolLoopAndStopMovement ()`

10.96.3.3 FindReference() [1/2] `void FindReference (`
 `unsigned char busaddress,`
 `char axes)`

10.96.3.4 FindReference() [2/2] `void FindReference (`
 `unsigned char busaddress,`
 `char axes,`
 `int timeout)`

Searches the reference position of the motor

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

10.96.3.5 GetAirpressure() `unsigned int GetAirpressure ()`

10.96.3.6 GetAirpressureLimit() `void GetAirpressureLimit (`
 `[System::Runtime::InteropServices::Out] unsigned int% lowerpressure,`
 `[System::Runtime::InteropServices::Out] unsigned int% upperpressure)`

10.96.3.7 GetAirValve() `unsigned int GetAirValve ()`

10.96.3.8 GetCurrentAirvalve() `unsigned int GetCurrentAirvalve ()`

10.96.3.9 GetCurrentAirvalveLimit() `void GetCurrentAirvalveLimit (`
 `[System::Runtime::InteropServices::Out] unsigned int% lowercurrent,`
 `[System::Runtime::InteropServices::Out] unsigned int% uppercurrent)`

10.96.3.10 GetCurrentPosition() `void GetCurrentPosition (`
 `unsigned char busaddress,`
 `char axes,`
 `[System::Runtime::InteropServices::Out] int% x,`
 `[System::Runtime::InteropServices::Out] int% y)`

Gets the current position of motors

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>x</i>	Current position of first axis if pattern in axes is set
<i>y</i>	Current position of second axis if pattern in axes is set

10.96.3.11 GetErrorAirpressure() `unsigned int GetErrorAirpressure ()`

10.96.3.12 GetErrorCurrentAirvalve() `unsigned int GetErrorCurrentAirvalve ()`

10.96.3.13 GetErrorVoltage12V() `unsigned int GetErrorVoltage12V ()`

10.96.3.14 GetErrorVoltage5V() `unsigned int GetErrorVoltage5V ()`

10.96.3.15 GetErrorVoltageAirvalve() `unsigned int GetErrorVoltageAirvalve ()`

10.96.3.16 GetErrorVoltageRs485A() unsigned int GetErrorVoltageRs485A ()

10.96.3.17 GetErrorVoltageRs485B() unsigned int GetErrorVoltageRs485B ()

10.96.3.18 GetErrorVoltageValves() unsigned int GetErrorVoltageValves ()

10.96.3.19 GetInMovement() bool GetInMovement ()

Low level command, gets the internal state "In Movement"

10.96.3.20 GetMinPressure() int GetMinPressure ()

10.96.3.21 GetMovementError() uint32_t GetMovementError ()

Low level command, gets the error of the last movement end

10.96.3.22 GetVoltage12V() unsigned int GetVoltage12V ()

10.96.3.23 GetVoltage12VLimit() void GetVoltage12VLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.24 GetVoltage5V() unsigned int GetVoltage5V ()

10.96.3.25 GetVoltage5VLimit() void GetVoltage5VLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.26 GetVoltageAirvalve() unsigned int GetVoltageAirvalve ()

10.96.3.27 GetVoltageAirvalveLimit() void GetVoltageAirvalveLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.28 GetVoltageRs485A() unsigned int GetVoltageRs485A ()

10.96.3.29 GetVoltageRs485ALimit() void GetVoltageRs485ALimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.30 GetVoltageRs485B() unsigned int GetVoltageRs485B ()

10.96.3.31 GetVoltageRs485BLimit() void GetVoltageRs485BLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.32 GetVoltageValves() unsigned int GetVoltageValves ()

10.96.3.33 GetVoltageValvesLimit() void GetVoltageValvesLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

10.96.3.34 MoveAbs() [1/2] void MoveAbs (
 unsigned char *busaddress*,
 char *axes*,
 int *x*,
 int *y*)

10.96.3.35 MoveAbs() [2/2] void MoveAbs (
 unsigned char *busaddress*,
 char *axes*,
 int *x*,
 int *y*,
 int *timeout*)

Moves the motor to the new absolute position

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>x</i>	Position of first axis, if pattern in axes is set
<i>y</i>	Position of second axis if pattern in axes is set
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

10.96.3.36 SetAirpressureLimit() `void SetAirpressureLimit (`
 `unsigned int lowerpressure,`
 `unsigned int upperpressure)`

10.96.3.37 SetAirValve() `void SetAirValve (`
 `unsigned int onoff)`

10.96.3.38 SetCurrentAirvalveLimit() `void SetCurrentAirvalveLimit (`
 `unsigned int lowercurrent,`
 `unsigned int uppercurrent)`

10.96.3.39 SetCurrentAndAir() [1/2] `void SetCurrentAndAir (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned short onoff)`

10.96.3.40 SetCurrentAndAir() [2/2] `void SetCurrentAndAir (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned short onoff,`
 `int timeout)`

10.96.3.41 SetInMovement() `void SetInMovement ()`

Low level command, sets the internal state to "In Movement"

10.96.3.42 SetMinPressure() void SetMinPressure (
int *pressure*)

10.96.3.43 SetVoltage12VLimit() void SetVoltage12VLimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.44 SetVoltage5VLimit() void SetVoltage5VLimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.45 SetVoltageAirvalveLimit() void SetVoltageAirvalveLimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.46 SetVoltageRs485ALimit() void SetVoltageRs485ALimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.47 SetVoltageRs485BLimit() void SetVoltageRs485BLimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.48 SetVoltageValvesLimit() void SetVoltageValvesLimit (
unsigned int *lowervoltage*,
unsigned int *uppervoltage*)

10.96.3.49 StopMovement() [1/2] void StopMovement (
unsigned char *busaddress*,
char *axes*)

10.96.3.50 StopMovement() [2/2] void StopMovement (
unsigned char *busaddress*,
char *axes*,
int *timeout*)

Stops the current movement

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

10.96.4 Member Data Documentation

10.96.4.1 Axes_I `const char Axes_I = 2 [static]`

Bit pattern for i axis for axes argument

10.96.4.2 Axes_X `const char Axes_X = 1 [static]`

Bit pattern for x axis for axes argument

10.96.4.3 Axes_Y `const char Axes_Y = 2 [static]`

Bit pattern for y axis for axes argument

10.96.4.4 Axes_Z `const char Axes_Z = 1 [static]`

Bit pattern for z axis for axes argument

10.96.4.5 Axis_I `const byte Axis_I = 1 [static]`

Axis number of i for axis argument

10.96.4.6 Axis_X `const byte Axis_X = 0 [static]`

Axis number of x for axis argument

10.96.4.7 Axis_Y `const byte Axis_Y = 1 [static]`

Axis number of y for axis argument

10.96.4.8 Axis_Z `const byte Axis_Z = 0 [static]`

Axis number of z for axis argument

10.96.4.9 McsBus_XY `const byte McsBus_XY = 1 [static]`

McsBus address for the xy-plane

10.96.4.10 McsBus_ZI `const byte McsBus_ZI = 2 [static]`

McsBus address for the z and i axes

10.96.4.11 RoboError_AnotherMaster `const uint32_t RoboError_AnotherMaster = ((0xA0110000L) | 0x0005) [static]`

10.96.4.12 RoboError_Base `const uint32_t RoboError_Base = (0xA0110000L) [static]`

10.96.4.13 RoboError_CannotEscapeEndSwitch `const uint32_t RoboError_CannotEscapeEndSwitch = ((0xA0110000L) | 0x0009) [static]`

10.96.4.14 RoboError_CommandAlreadyInProgress `const uint32_t RoboError_CommandAlreadyInProgress = ((0xA0110000L) | 0x000A) [static]`

10.96.4.15 RoboError_CommandNotPossible `const uint32_t RoboError_CommandNotPossible = ((0xA0110000L) | 0x0014) [static]`

10.96.4.16 RoboError_CommunicationTimeout `const uint32_t RoboError_CommunicationTimeout = ((0xA0110000L) | 0x0004) [static]`

10.96.4.17 RoboError_DacqNotReady `const uint32_t RoboError_DacqNotReady = ((0xA0110000L) | 0x0015) [static]`

10.96.4.18 RoboError_DLLMovementTimeout `const uint32_t RoboError_DLLMovementTimeout = ((0xA0110000L) | 0x1001) [static]`

10.96.4.19 RoboError_FindReferenceMethod `const uint32_t RoboError_FindReferenceMethod = ((0xA0110000L) | 0x0006) [static]`

10.96.4.20 RoboError_GilsonCommandPending `const uint32_t RoboError_GilsonCommandPending = ((0xA0110000L) | 0x0011) [static]`

10.96.4.21 RoboError_GilsonTimeout `const uint32_t RoboError_GilsonTimeout = ((0xA0110000L) | 0x000F) [static]`

10.96.4.22 RoboError_GilsonWrondID `const uint32_t RoboError_GilsonWrondID = ((0xA0110000L) | 0x0010) [static]`

10.96.4.23 RoboError_McsBus_UnknownCommand `const uint32_t RoboError_McsBus_UnknownCommand = ((0xA0110000L) | 0x003F) [static]`

10.96.4.24 RoboError_NoEndSwitch `const uint32_t RoboError_NoEndSwitch = ((0xA0110000L) | 0x0008) [static]`

10.96.4.25 RoboError_NoMoreData `const uint32_t RoboError_NoMoreData = ((0xA0110000L) | 0x0016) [static]`

10.96.4.26 RoboError_NoReference `const uint32_t RoboError_NoReference = ((0xA0110000L) | 0x000B) [static]`

10.96.4.27 RoboError_NoSpeedOrAcceleration `const uint32_t RoboError_NoSpeedOrAcceleration = ((0xA0110000L) | 0x0007) [static]`

10.96.4.28 RoboError_OverPressure `const uint32_t RoboError_OverPressure = ((0xA0110000L) | 0x000C) [static]`

10.96.4.29 RoboError_ParameterNotAllowed `const uint32_t RoboError_ParameterNotAllowed = ((0xA0110000L) | 0x0012) [static]`

10.96.4.30 RoboError_PeristalticTimeout `const uint32_t RoboError_PeristalticTimeout = ((0xA0110000L) | 0x000E) [static]`

10.96.4.31 RoboError_Phase0OutOfRange `const uint32_t RoboError_Phase0OutOfRange = ((0xA0110000L) | 0x000D) [static]`

10.96.4.32 RoboError_PollLoopCanceled `const uint32_t RoboError_PollLoopCanceled = ((0xA0110000L) | 0x1002) [static]`

10.96.4.33 RoboError_PollLoopCanceledAndStopMovement `const uint32_t RoboError_PollLoopCanceledAndStopMovement = ((0xA0110000L) | 0x1003) [static]`

10.96.4.34 RoboError_Pressure `const uint32_t RoboError_Pressure = ((0xA0110000L) | 0x0002) [static]`

10.96.4.35 RoboError_RangeExceeded `const uint32_t RoboError_RangeExceeded = ((0xA0110000L) | 0x0003) [static]`

10.96.4.36 RoboError_StateChangeNotPossible `const uint32_t RoboError_StateChangeNotPossible = ((0xA0110000L) | 0x0013) [static]`

10.96.4.37 RoboError_Timeout `const uint32_t RoboError_Timeout = ((0xA0110000L) | 0x0001) [static]`

10.96.4.38 RoboError_UnknownCommand `const uint32_t RoboError_UnknownCommand = ((0xA0110000L)) [static]`

10.96.5 Property Documentation

10.96.5.1 McsBus `CMcsBusNet^ McsBus [get]`

10.96.5.2 McsBus_MotorControl `CMcsBus_MotorControlNet^ McsBus_MotorControl [get]`


10.96.5.3 RoboMainLowLevelCommand `RoboMainLowLevelCommands^ RoboMainLowLevelCommand [get]`

10.96.6 Event Documentation

10.96.6.1 RoboStatusEvent `RoboStatusEventDelegate^ RoboStatusEvent`

10.97 CRoboFluidDeviceNet Class Reference

Inheritance diagram for CRoboFluidDeviceNet:

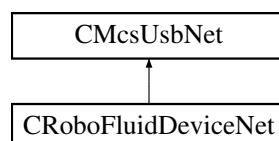

Public Member Functions

- [CRoboFluidDeviceNet](#) (void)
- [~CRoboFluidDeviceNet](#) (void)
- void [SetValve](#) (int value)
Open or Close valves.
- void [SetSingleValve](#) (int valve, bool onoff)
Opens or Closes a valve.
- int [GetValve](#) ()
Query the state of the values.
- bool [GetSingleValve](#) (int valve)
Query the state of a valve.
- void [CloseAllValves](#) ()
- void [PumpOn](#) (int index, short speed)
- void [SetPumpSpeed](#) (int index, short speed)
- void [PumpOff](#) (int index)
- short [GetPumpSpeed](#) (int index)
- bool [IsPumpMotorOn](#) (int index)

Protected Attributes

- CRoboFluidDevice * [m_pRoboFluidDevice](#)
- [CMcsBus_MotorControlNet](#) ^ [m_pMcsBus_MotorControlNet](#)

Properties

- [CMcsBus_MotorControlNet](#) ^ [McsBus_MotorControl](#) [get]

Additional Inherited Members

10.97.1 Constructor & Destructor Documentation

10.97.1.1 CRoboFluidDeviceNet() [CRoboFluidDeviceNet](#) (
void)

10.97.1.2 ~CRoboFluidDeviceNet() [~CRoboFluidDeviceNet](#) (
void)

10.97.2 Member Function Documentation

10.97.2.1 CloseAllValves() `void CloseAllValves ()`

10.97.2.2 GetPumpSpeed() `short GetPumpSpeed (`
`int index)`

10.97.2.3 GetSingleValve() `bool GetSingleValve (`
`int valve)`

Query the state of a valve.

Parameters

<i>valve</i>	number of valve /*!
--------------	---------------------

Returns

state of the valve

10.97.2.4 GetValve() `int GetValve ()`

Query the state of the values.

Returns

the current state of the valves as a bit pattern.

10.97.2.5 IsPumpMotorOn() `bool IsPumpMotorOn (`
`int index)`

10.97.2.6 PumpOff() `void PumpOff (`
`int index)`

10.97.2.7 PumpOn() `void PumpOn (`
`int index,`
`short speed)`

10.97.2.8 SetPumpSpeed() `void SetPumpSpeed (`
 `int index,`
 `short speed)`

10.97.2.9 SetSingleValve() `void SetSingleValve (`
 `int valve,`
 `bool onoff)`

Opens or Closes a valve.

Parameters

<i>valve</i>	number of valve to be changed /*!
<i>onoff</i>	open or close the valve

10.97.2.10 SetValve() `void SetValve (`
 `int value)`

Open or Close valves.

Parameters

<i>value</i>	bit pattern of valves which should be open.
--------------	---

10.97.3 Member Data Documentation

10.97.3.1 m_pMcsBus_MotorControlNet `CMcsBus_MotorControlNet ^ m_pMcsBus_MotorControlNet`
 [protected]

10.97.3.2 m_pRoboFluidDevice `CRoboFluidDevice* m_pRoboFluidDevice` [protected]


10.97.4 Property Documentation

10.97.4.1 McsBus_MotorControl `CMcsBus_MotorControlNet^ McsBus_MotorControl` [get]

10.98 CRobolInjectDeviceNet Class Reference

[CRobolInjectDeviceNet](#) is the to control the MCS RobolInject device

Inheritance diagram for CRobolInjectDeviceNet:

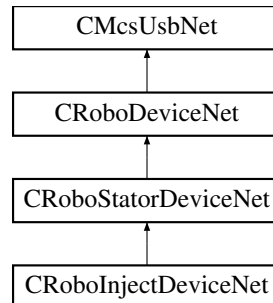

Public Member Functions

- [CRobolInjectDeviceNet](#) (void)

Additional Inherited Members

10.98.1 Detailed Description

[CRobolInjectDeviceNet](#) is the to control the MCS RobolInject device


10.98.2 Constructor & Destructor Documentation

10.98.2.1 CRobolInjectDeviceNet() [CRobolInjectDeviceNet](#) (
void)

10.99 CRoboocyte2DeviceNet Class Reference

[CRoboocyte2DeviceNet](#) is the class to control the MCS Roboocyte2 device

Inheritance diagram for CRoboocyte2DeviceNet:

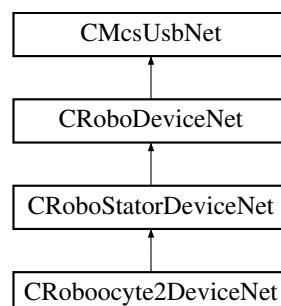

Public Member Functions

- [CRobocyte2DeviceNet](#) (void)
- void [SetAxisLED](#) (bool onoff)
- bool [GetAxisLED](#) ()
- [CRoboDacqNet](#) ^ [GetRoboDacq](#) ()
- [CRoboFluidDeviceNet](#) ^ [GetRoboFluidDevice](#) ()
- [CGilsonDeviceNet](#) ^ [GetGilsonDevice](#) ()
- [CMcsBus_ExtensionNet](#) ^ [GetMcsBus_Extension](#) ()

Additional Inherited Members

10.99.1 Detailed Description

[CRobocyte2DeviceNet](#) is the class to control the MCS Robocyte2 device

10.99.2 Constructor & Destructor Documentation

10.99.2.1 CRobocyte2DeviceNet() [CRobocyte2DeviceNet](#) (
void)

10.99.3 Member Function Documentation

10.99.3.1 GetAxisLED() bool [GetAxisLED](#) ()

10.99.3.2 GetGilsonDevice() [CGilsonDeviceNet](#) ^ [GetGilsonDevice](#) ()

10.99.3.3 GetMcsBus_Extension() [CMcsBus_ExtensionNet](#) ^ [GetMcsBus_Extension](#) ()


10.99.3.4 GetRoboDacq() [CRoboDacqNet](#) ^ [GetRoboDacq](#) ()

10.99.3.5 GetRoboFluidDevice() [CRoboFluidDeviceNet](#) ^ GetRoboFluidDevice ()

10.99.3.6 SetAxisLED() void SetAxisLED (
 bool onoff)

10.100 CRoboStatorDeviceNet Class Reference

Inheritance diagram for CRoboStatorDeviceNet:

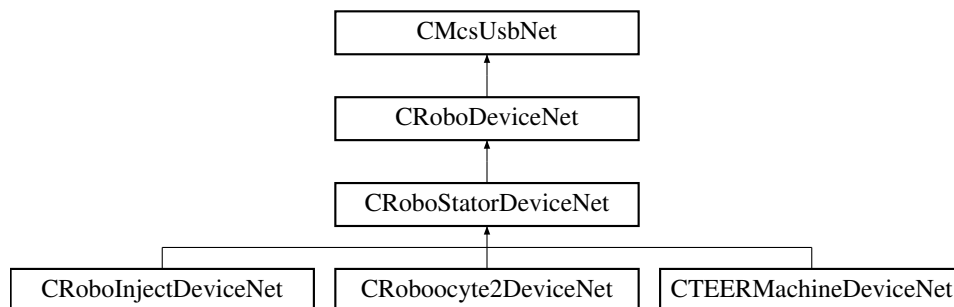

Classes

- class [RoboMainStatorLowLevelCommands](#)

Public Member Functions

- [CRoboStatorDeviceNet](#) (void)
- void [FindReferenceXY](#) ()
- void [FindReferenceXY](#) (int timeout)
- void [FindReferenceZ](#) ()
- void [FindReferenceZ](#) (int timeout)
- void [FindReferenceI](#) ()
- void [FindReferenceI](#) (int timeout)
- unsigned char [HasRefXY](#) ()
- unsigned char [HasRefZ](#) ()
- unsigned char [HasRefI](#) ()
- void [MoveAbsXY](#) (int x, int y)
- void [MoveAbsXY](#) (int x, int y, int timeout)
- void [MoveAbsZ](#) (int z)
- void [MoveAbsZ](#) (int z, int timeout)
- void [MoveAbsI](#) (int i)
- void [MoveAbsI](#) (int i, int timeout)
- void [StopMovementXY](#) ()
- void [StopMovementXY](#) (int timeout)
- void [StopMovementZ](#) ()
- void [StopMovementZ](#) (int timeout)
- void [StopMovementI](#) ()
- void [StopMovementI](#) (int timeout)
- void [SetCurrentAndAirXY](#) (unsigned short onoff)

- void [SetCurrentAndAirXY](#) (unsigned short onoff, int timeout)
- void [GetCurrentPositionXY](#) ([System::Runtime::InteropServices::Out]int% x, [System::Runtime::InteropServices::Out]int% y)
- void [GetCurrentPositionZ](#) ([System::Runtime::InteropServices::Out]int% z)
- void [GetCurrentPositionI](#) ([System::Runtime::InteropServices::Out]int% i)
- void [SetVelocityXY](#) (int v)
- void [SetVelocityZ](#) (int v)
- void [SetVelocityI](#) (int v)
- void [SetSpeedXY](#) (int v)
- void [SetSpeedZ](#) (int v)
- void [SetSpeedI](#) (int v)
- void [SetSpeedNativeXY](#) (int v)
- void [SetSpeedNativeZ](#) (int v)
- void [SetSpeedNativeI](#) (int v)
- void [SetAccelerationXY](#) (int a)
- void [SetAccelerationZ](#) (int a)
- void [SetAccelerationI](#) (int a)
- void [SetAccelerationNativeXY](#) (int a)
- void [SetAccelerationNativeZ](#) (int a)
- void [SetAccelerationNativeI](#) (int a)

Properties

- [RoboMainStatorLowLevelCommands](#)[^] [RoboMainStatorLowLevelCommand](#) [get]

Additional Inherited Members

10.100.1 Constructor & Destructor Documentation

10.100.1.1 CRoboStatorDeviceNet() [CRoboStatorDeviceNet](#) (
 void)

10.100.2 Member Function Documentation

10.100.2.1 FindReferenceI() [1/2] void FindReferenceI ()

10.100.2.2 FindReferenceI() [2/2] void FindReferenceI (
 int timeout)

10.100.2.3 FindReferenceXY() [1/2] void FindReferenceXY ()

10.100.2.4 FindReferenceXY() [2/2] void FindReferenceXY (
int *timeout*)

10.100.2.5 FindReferenceZ() [1/2] void FindReferenceZ ()

10.100.2.6 FindReferenceZ() [2/2] void FindReferenceZ (
int *timeout*)

10.100.2.7 GetCurrentPositionI() void GetCurrentPositionI (
[System::Runtime::InteropServices::Out] int% *i*)

10.100.2.8 GetCurrentPositionXY() void GetCurrentPositionXY (
[System::Runtime::InteropServices::Out] int% *x*,
[System::Runtime::InteropServices::Out] int% *y*)

10.100.2.9 GetCurrentPositionZ() void GetCurrentPositionZ (
[System::Runtime::InteropServices::Out] int% *z*)

10.100.2.10 HasRefI() unsigned char HasRefI ()

10.100.2.11 HasRefXY() unsigned char HasRefXY ()

10.100.2.12 HasRefZ() unsigned char HasRefZ ()

10.100.2.13 MoveAbsI() [1/2] void MoveAbsI (
int *i*)

10.100.2.14 MoveAbsI() [2/2] void MoveAbsI (
int *i*,
int *timeout*)

10.100.2.15 MoveAbsXY() [1/2] void MoveAbsXY (
int *x*,
int *y*)

10.100.2.16 MoveAbsXY() [2/2] void MoveAbsXY (
int *x*,
int *y*,
int *timeout*)

10.100.2.17 MoveAbsZ() [1/2] void MoveAbsZ (
int *z*)

10.100.2.18 MoveAbsZ() [2/2] void MoveAbsZ (
int *z*,
int *timeout*)

10.100.2.19 SetAccelerationI() void SetAccelerationI (
int *a*)

10.100.2.20 SetAccelerationNativeI() void SetAccelerationNativeI (
int *a*)

10.100.2.21 SetAccelerationNativeXY() void SetAccelerationNativeXY (
int *a*)

10.100.2.22 SetAccelerationNativeZ() `void SetAccelerationNativeZ (`
`int a)`

10.100.2.23 SetAccelerationXY() `void SetAccelerationXY (`
`int a)`

10.100.2.24 SetAccelerationZ() `void SetAccelerationZ (`
`int a)`

10.100.2.25 SetCurrentAndAirXY() [1/2] `void SetCurrentAndAirXY (`
`unsigned short onoff)`

10.100.2.26 SetCurrentAndAirXY() [2/2] `void SetCurrentAndAirXY (`
`unsigned short onoff,`
`int timeout)`

10.100.2.27 SetSpeedI() `void SetSpeedI (`
`int v)`

10.100.2.28 SetSpeedNativeI() `void SetSpeedNativeI (`
`int v)`

10.100.2.29 SetSpeedNativeXY() `void SetSpeedNativeXY (`
`int v)`

10.100.2.30 SetSpeedNativeZ() `void SetSpeedNativeZ (`
`int v)`

10.100.2.31 SetSpeedXY() void SetSpeedXY (
int v)

10.100.2.32 SetSpeedZ() void SetSpeedZ (
int v)

10.100.2.33 SetVelocityI() void SetVelocityI (
int v)

10.100.2.34 SetVelocityXY() void SetVelocityXY (
int v)

10.100.2.35 SetVelocityZ() void SetVelocityZ (
int v)

10.100.2.36 StopMovementI() [1/2] void StopMovementI ()

10.100.2.37 StopMovementI() [2/2] void StopMovementI (
int timeout)

10.100.2.38 StopMovementXY() [1/2] void StopMovementXY ()

10.100.2.39 StopMovementXY() [2/2] void StopMovementXY (
int timeout)

10.100.2.40 StopMovementZ() [1/2] void StopMovementZ ()


10.100.2.41 StopMovementZ() [2/2] void StopMovementZ (
int timeout)

10.100.3 Property Documentation

10.100.3.1 RoboMainStatorLowLevelCommand [RoboMainStatorLowLevelCommands](#)^ RoboMainStator↔
LowLevelCommand [get]

10.101 CSafeISDeviceNet Class Reference

Inheritance diagram for CSafeISDeviceNet:

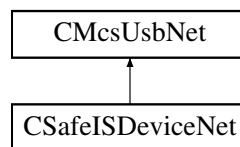

Public Member Functions

- [CSafeISDeviceNet](#) (void)
Initializes a new instance of the [CSafeISDeviceNet](#) class.
- [~CSafeISDeviceNet](#) (void)
Releases unmanaged resources and performs other cleanup operations before the [CSafeISDeviceNet](#) is reclaimed by garbage collection.
- void [SetSwitches](#) (unsigned short switches)
Sets the switches for all electrodes on the device. Do not use during measurement
- void [SetAdcChannels](#) (unsigned char channels)
Sets the ADC channels you want to be sampled
- void [SetAdcSamplePos](#) (array< unsigned short >^ positions)
Sets the sample position of the ADC.
- void [SetDacMode](#) (unsigned char mode)
Sets the DAC mode.
- void [SetDacPulseform](#) (array< short >^ pulseform)
Sets the DAC pulseform.
- void [SetDacPeriode](#) (unsigned int periode)
Sets the DAC periode.

Properties

- [CRoboDeviceNet](#)^ [RoboDevice](#) [get]
Gets the [CRoboDeviceNet](#). Use this to control the syringe.
- [CFluidControlDeviceNet](#)^ [FluidControlDevice](#) [get]
Gets the [CFluidControlDeviceNet](#). Use this to control the valves. Only SetSingleValve is implemented for [CSafeISDeviceNet](#).
- [CMcsUsbDacqNet](#)^ [DacqDevice](#) [get]
Gets the [CMcsUsbDacqNet](#). Use this to control the data aquisition.

Additional Inherited Members

10.101.1 Detailed Description

10.101.2 Constructor & Destructor Documentation

10.101.2.1 CSafeISDeviceNet() CSafeISDeviceNet (void)

Initializes a new instance of the [CSafeISDeviceNet](#) class.

10.101.2.2 ~CSafeISDeviceNet() ~CSafeISDeviceNet (void)

Releases unmanaged resources and performs other cleanup operations before the [CSafeISDeviceNet](#) is reclaimed by garbage collection.

10.101.3 Member Function Documentation

10.101.3.1 SetAdcChannels() void SetAdcChannels (unsigned char *channels*)

Sets the ADC channels you want to be sampled

Parameters

<i>channels</i>	The bitmap of the 8 channels. Set bit to 1 for the channels you want measure
-----------------	--

10.101.3.2 SetAdcSamplePos() void SetAdcSamplePos (array< unsigned short >^ *positions*)

Sets the sample position of the ADC.

Parameters

<i>positions</i>	The positions in units of 0.1μs.
------------------	----------------------------------

10.101.3.3 SetDacMode() `void SetDacMode (`
`unsigned char mode)`

Sets the DAC mode.

Parameters

<i>mode</i>	The mode: 0 = Impedance ; 1 = Amperometry
-------------	---

10.101.3.4 SetDacPeriode() `void SetDacPeriode (`
`unsigned int periode)`

Sets the DAC periode.

Parameters

<i>periode</i>	The periode in units of 10µs.
----------------	-------------------------------

10.101.3.5 SetDacPulseform() `void SetDacPulseform (`
`array< short >^ pulseform)`

Sets the DAC pulseform.

Parameters

<i>pulseform</i>	The pulseform.
------------------	----------------

10.101.3.6 SetSwitches() `void SetSwitches (`
`unsigned short switches)`

Sets the switches for all electrodes on the device. Do not use during measurement

Parameters

<i>switches</i>	The switches: See Schematics for the meaning
-----------------	--

10.101.4 Property Documentation

10.101.4.1 DacqDevice `CMcsUsbDacqNet^ DacqDevice [get]`

Gets the `CMcsUsbDacqNet`. Use this to control the data acquisition.

10.101.4.2 FluidControlDevice `CFluidControlDeviceNet^ FluidControlDevice [get]`


Gets the `CFluidControlDeviceNet`. Use this to control the valves. Only `SetSingleValve` is implemented for `CSafeISDeviceNet`.

10.101.4.3 RoboDevice `CRoboDeviceNet^ RoboDevice [get]`

Gets the `CRoboDeviceNet`. Use this to control the syringe.

10.102 CSCUDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for `CSCUDacqGroupChannelSelectionNet`:

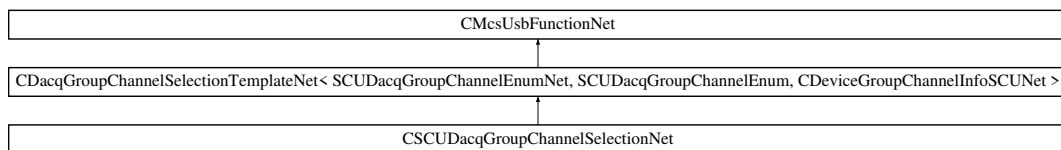**Public Member Functions**


- `CSCUDacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)`

Additional Inherited Members**10.102.1 Constructor & Destructor Documentation****10.102.1.1 CSCUDacqGroupChannelSelectionNet()** `CSCUDacqGroupChannelSelectionNet (
 CMcsUsbNet^ mcsusb)`

10.103 CSCUFunctionNet Class Reference

[CSCUFunctionNet](#) is the class to control the SCU device

Inheritance diagram for CSCUFunctionNet:

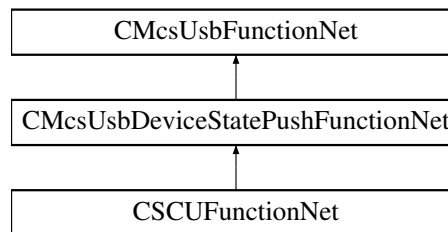

Public Member Functions

- delegate void [OnGetAvailableHeadstages](#) (uint32_t AvailableHeadstages)
- delegate void [OnIsHeadstageAvailable](#) (uint32_t Headstage, bool available)
- [CSCUFunctionNet](#) (CMcsUsbNet[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pSCUFunctionPointerContainer)
 - Initializes a new instance of the [CSCUFunctionNet](#) class.*
- [CSCUFunctionNet](#) (CMcsUsbNet[^] mcsusb)
- virtual [~CSCUFunctionNet](#) ()
- [!CSCUFunctionNet](#) ()
- uint32_t [GetAvailableHeadstages](#) ()
 - Gets a bitmap of available headstages.*
- bool [IsInDacqLegacyMode](#) ()
 - Is the SCU in legacy mode*
- void [SetDacqLegacyMode](#) (bool enable)
 - Enable the SCU legacy mode*
- uint32_t [GetMaxStimulusChannelsPerHeadstage](#) ()
 - Gets the maximal number of stimulation channels a headstage can have.*
- uint32_t [GetMaxNumberOfHeadstages](#) ()
 - Gets the maximal number of headstages.*
- SCU_HeadstageIdEnumNet [GetHeadstageID](#) (uint32_t Headstage)
 - Gets the headstage fpga ID.*
- bool [IsHeadstageAvailable](#) (uint32_t Headstage)
 - Checks whether the given headstage is available.*
- void [PowerHS](#) (uint32_t Headstage, bool power)
 - Power the HS*
- bool [IsHSPowered](#) (uint32_t Headstage)
 - Is the HS powered*
- bool [HasHSPowerSwitch](#) ()
 - Has SCU HS power switch*
- String[^] [GetHeadstageSerialNumber](#) (uint32_t Headstage)
 - Gets the serial number of a given headstage.*
- uint32_t [GetHeadstageNumberOfAnalogChannels](#) (uint32_t Headstage)
 - Gets the number of analog channels for a given headstage.*
- uint32_t [GetHeadstageNumberOfStimulationChannels](#) (uint32_t Headstage)
 - Gets the number of stimulation channels for a given headstage.*
- uint32_t [GetHeadstageGainInPer mille](#) (uint32_t Headstage)

- Gets the gain factor in permille for a given headstage.*
- uint32_t [GetHeadstageAdcRangeInMicroVolt](#) (uint32_t Headstage)
 - Gets the ADC Range in uV for a given headstage.*
- uint32_t [GetHeadstageAdcBits](#) (uint32_t Headstage)
 - Gets the Number of ADC bits for a given headstage.*
- uint32_t [GetHeadstageDacVoltageRangeInMilliVolt](#) (uint32_t Headstage)
 - Gets the DAC Voltage Range in mV for a given headstage.*
- uint32_t [GetHeadstageDacVoltageResolutionInMicroVolt](#) (uint32_t Headstage)
 - Gets the DAC Voltage Resolution in uV for a given headstage.*
- uint32_t [GetHeadstageDacCurrentRangeInMicroAmpere](#) (uint32_t Headstage)
 - Gets the DAC Current Range in uA for a given headstage.*
- uint32_t [GetHeadstageDacCurrentResolutionInNanoAmpere](#) (uint32_t Headstage)
 - Gets the DAC Current Resolution in nA for a given headstage.*
- uint32_t [GetHeadstageDacBits](#) (uint32_t Headstage)
 - Gets the Number of DAC bits for a given headstage.*
- uint32_t [GetHeadstageSamplerate](#) (uint32_t Headstage)
 - Gets the Samplerate of a given headstage.*
- bool [HasGalvanicIsolation](#) ()
 - Has galvanic isolated hardware*
- bool [HasAnalogOut](#) ()
 - Has AnalogOut hardware*
- void [EnableAnalogOut](#) (bool enable)
 - Enables AnalogOut globally*
- bool [IsAnalogOutEnabled](#) ()
 - Is AnalogOut enabled*
- void [SetAnalogOutDACRange](#) (AnalogOut_DAC_Range_EnumNet range)
 - Sets the analog out DAC range*
- AnalogOut_DAC_Range_EnumNet [GetAnalogOutDACRange](#) ()
 - Gets the analog out DAC range*
- void [SetAnalogOutADCRange](#) (uint32_t range)
 - Sets the analog out ADC range*
- uint32_t [GetAnalogOutADCRange](#) ()
 - Gets the analog out ADC range*
- void [AutomaticAnalogOut](#) (bool automatic)
 - Sets automatic source channel selection*
- bool [IsAutomaticAnalogOut](#) ()
 - Is Automatic source channel selection selected*
- void [SetAnalogOutChannels](#) (uint32_t out_channel, uint32_t source_channel)
 - Set the source channel number for a certain output channel*
- uint32_t [GetAnalogOutChannels](#) (uint32_t out_channel)
 - Get the connected source channel number for a certain output channel*
- void [SetReferenceElectrodeSwitchState](#) (uint32_t Headstage, ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)
 - Sets the position of the switch for the reference electrode*
- ReferenceElectrodeSwitchPositionEnumNet [GetReferenceElectrodeSwitchState](#) (uint32_t Headstage)
 - Gets the position of the switch for the reference electrode*
- void [SetReferenceElectrodeMode](#) (uint32_t Headstage, ReferenceElectrodeModeEnumNet NewValue)
 - Sets the mode for the reference electrode*
- ReferenceElectrodeModeEnumNet [GetReferenceElectrodeMode](#) (uint32_t Headstage)
 - Gets the mode for the reference electrode*
- [CFilterPropertyNet](#) ^ [GetFilterProperty](#) (SCUDacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)

Gets the filter property

- array< [CFilterPropertyNet](#)> ^ [GetFilterProperties](#) (SCUDacqGroupChannelEnumNet GroupID, int filter← Configurations_Length)

Gets multiple filter properties

Events

- [OnGetAvailableHeadstages](#) ^ [GetAvailableHeadstagesEvent](#) [add, remove, raise]
Event fires when the bitmap of available headstages has changed
- [OnIsHeadstageAvailable](#) ^ [IsHeadstageAvailableEvent](#) [add, remove, raise]
Event fires when 'true' if the headstage is connected for the headstage to query has changed

Additional Inherited Members

10.103.1 Detailed Description

[CSCUFunctionNet](#) is the class to control the SCU device

10.103.2 Constructor & Destructor Documentation

10.103.2.1 [CSCUFunctionNet\(\)](#) [1/2] [CSCUFunctionNet](#) (
 [CMcsUsbNet](#) ^ *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#) ^ *pSCUFunctionPointerContainer*)

Initializes a new instance of the [CSCUFunctionNet](#) class.

10.103.2.2 [CSCUFunctionNet\(\)](#) [2/2] [CSCUFunctionNet](#) (
 [CMcsUsbNet](#) ^ *mcsusb*)

10.103.2.3 [~CSCUFunctionNet\(\)](#) virtual [~CSCUFunctionNet](#) () [virtual]

10.103.2.4 ["!CSCUFunctionNet\(\)](#) [!CSCUFunctionNet](#) ()

10.103.3 Member Function Documentation

10.103.3.1 [AutomaticAnalogOut\(\)](#) void [AutomaticAnalogOut](#) (
 bool *automatic*)

Sets automatic source channel selection

Parameters

<i>automatic</i>	Automatic
------------------	-----------

10.103.3.2 EnableAnalogOut() `void EnableAnalogOut (`
 `bool enable)`

Enables AnalogOut globally

Parameters

<i>enable</i>	Enable
---------------	--------

10.103.3.3 GetAnalogOutADCRange() `uint32_t GetAnalogOutADCRange ()`

Gets the analog out ADC range

Returns

Range

10.103.3.4 GetAnalogOutChannels() `uint32_t GetAnalogOutChannels (`
 `uint32_t out_channel)`

Get the connected source channel number for a certain output channel

Parameters

<i>out_channel</i>	Output channel number
--------------------	-----------------------

Returns

Source channel number

10.103.3.5 GetAnalogOutDACRange() `AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ()`

Gets the analog out DAC range

Returns

Range

10.103.3.6 GetAvailableHeadstages() `uint32_t GetAvailableHeadstages ()`

Gets a bitmap of available headstages.

Returns

The bitmap of available headstages.

10.103.3.7 GetFilterProperties() `array<CFilterPropertyNet^> ^ GetFilterProperties (SCUDacqGroupChannelEnumNet GroupID, int filterConfigurations_Length)`

Gets multiple filter properties

Parameters

<i>GroupID</i>	The group ID
<i>filterConfigurations_Length</i>	The maximal length of filterConfigurations.

Returns

array of filter properties

10.103.3.8 GetFilterProperty() `CFilterPropertyNet ^ GetFilterProperty (SCUDacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`

Gets the filter property

Parameters

<i>GroupID</i>	The group ID
<i>FilterNumber</i>	The filter number

Returns

The filter property

10.103.3.9 GetHeadstageAdcBits() `uint32_t GetHeadstageAdcBits (uint32_t Headstage)`

Gets the Number of ADC bits for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of bits the ADC has for the given headstage.

10.103.3.10 GetHeadstageAdcRangeInMicroVolt() `uint32_t GetHeadstageAdcRangeInMicroVolt (uint32_t Headstage)`

Gets the ADC Range in uV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The ADC Range in uV for the given headstage.

10.103.3.11 GetHeadstageDacBits() `uint32_t GetHeadstageDacBits (uint32_t Headstage)`

Gets the Number of DAC bits for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of bits the DAC has for the given headstage.

10.103.3.12 GetHeadstageDacCurrentRangeInMicroAmpere() `uint32_t GetHeadstageDacCurrentRangeInMicroAmpere (uint32_t Headstage)`

Gets the DAC Current Range in uA for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Current Range in uA for the given headstage.

10.103.3.13 GetHeadstageDacCurrentResolutionInNanoAmpere() uint32_t GetHeadstageDacCurrentResolutionInNanoAmpere (
 uint32_t Headstage)

Gets the DAC Current Resolution in nA for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Current Resolution in nA for the given headstage.

10.103.3.14 GetHeadstageDacVoltageRangeInMilliVolt() uint32_t GetHeadstageDacVoltageRangeInMilliVolt (
 uint32_t Headstage)

Gets the DAC Voltage Range in mV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Voltage Range in mV for the given headstage.

10.103.3.15 GetHeadstageDacVoltageResolutionInMicroVolt() uint32_t GetHeadstageDacVoltageResolutionInMicroVolt (
 uint32_t Headstage)

Gets the DAC Voltage Resolution in uV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Voltage Resolution in uV for the given headstage.

10.103.3.16 GetHeadstageGainInPer mille() `uint32_t GetHeadstageGainInPer mille (uint32_t Headstage)`

Gets the gain factor in permille for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The gain factor in permille for the given headstage.

10.103.3.17 GetHeadstageID() `SCU_HeadstageIdEnumNet GetHeadstageID (uint32_t Headstage)`

Gets the headstage fpga ID.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The headstage fpga ID.

10.103.3.18 GetHeadstageNumberOfAnalogChannels() `uint32_t GetHeadstageNumberOfAnalogChannels (uint32_t Headstage)`

Gets the number of analog channels for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of analog channels the headstage has.

10.103.3.19 GetHeadstageNumberOfStimulationChannels() `uint32_t GetHeadstageNumberOfStimulationChannels (uint32_t Headstage)`

Gets the number of stimulation channels for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of stimulation channels the headstage has.

10.103.3.20 GetHeadstageSamplerate() `uint32_t GetHeadstageSamplerate (uint32_t Headstage)`

Gets the Samplerate of a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The samplerate in Hz for the given headstage.

10.103.3.21 GetHeadstageSerialNumber() `String ^ GetHeadstageSerialNumber (uint32_t Headstage)`

Gets the serial number of a given headstage.

10.103.3.25 GetReferenceElectrodeSwitchState() `ReferenceElectrodeSwitchPositionEnumNet Get↔
ReferenceElectrodeSwitchState (
 uint32_t Headstage)`

Gets the position of the switch for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The switch position

10.103.3.26 HasAnalogOut() `bool HasAnalogOut ()`

Has AnalogOut hardware

Returns

Enabled

10.103.3.27 HasGalvanicIsolation() `bool HasGalvanicIsolation ()`

Has galvanic isolated hardware

Returns

Enabled

10.103.3.28 HasHSPowerSwitch() `bool HasHSPowerSwitch ()`

Has SCU HS power switch

Returns

Has Switch

10.103.3.29 IsAnalogOutEnabled() `bool IsAnalogOutEnabled ()`

Is AnalogOut enabled

Returns

Enabled

10.103.3.30 IsAutomaticAnalogOut() `bool IsAutomaticAnalogOut ()`

Is Automatic source channel selection selected

Returns

Automatic

10.103.3.31 IsHeadstageAvailable() `bool IsHeadstageAvailable (
uint32_t Headstage)`

Checks whether the given headstage is available.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

'true' if the headstage is connected.

10.103.3.32 IsHSPowered() `bool IsHSPowered (
uint32_t Headstage)`

Is the HS powered

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

'true' if the headstage is powered.

10.103.3.33 IsInDacqLegacyMode() `bool IsInDacqLegacyMode ()`

Is the SCU in legacy mode

Returns

Is Enabled

10.103.3.34 OnGetAvailableHeadstages() `delegate void OnGetAvailableHeadstages (
uint32_t AvailableHeadstages)`**10.103.3.35 OnIsHeadstageAvailable()** `delegate void OnIsHeadstageAvailable (
uint32_t Headstage,
bool available)`**10.103.3.36 PowerHS()** `void PowerHS (
uint32_t Headstage,
bool power)`

Power the HS

Parameters

<i>Headstage</i>	The headstage to query.
<i>power</i>	'true' if the headstage is powered.

10.103.3.37 SetAnalogOutADCRange() `void SetAnalogOutADCRange (
uint32_t range)`

Sets the analog out ADC range

Parameters

<i>range</i>	Range
--------------	-------

10.103.3.38 SetAnalogOutChannels() `void SetAnalogOutChannels (
uint32_t out_channel,
uint32_t source_channel)`

Set the source channel number for a certain output channel

Parameters

<i>out_channel</i>	Output channel number
<i>source_channel</i>	Source channel number

10.103.3.39 SetAnalogOutDACRange() `void SetAnalogOutDACRange (AnalogOut_DAC_Range_EnumNet range)`

Sets the analog out DAC range

Parameters

<i>range</i>	Range
--------------	-------

10.103.3.40 SetDacqLegacyMode() `void SetDacqLegacyMode (bool enable)`

Enable the SCU legacy mode

Parameters

<i>enable</i>	Enable
---------------	--------

10.103.3.41 SetReferenceElectrodeMode() `void SetReferenceElectrodeMode (uint32_t Headstage, ReferenceElectrodeModeEnumNet NewValue)`

Sets the mode for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
<i>NewValue</i>	The mode

10.103.3.42 SetReferenceElectrodeSwitchState() `void SetReferenceElectrodeSwitchState (uint32_t Headstage, ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)`

Sets the position of the switch for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
<i>NewSwitchPos</i>	The switch position

10.103.4 Event Documentation

10.103.4.1 GetAvailableHeadstagesEvent [OnGetAvailableHeadstages](#)[^] [GetAvailableHeadstagesEvent](#) [add], [remove], [raise]


Event fires when the bitmap of available headstages has changed

10.103.4.2 IsHeadstageAvailableEvent [OnIsHeadstageAvailable](#)[^] [IsHeadstageAvailableEvent](#) [add], [remove], [raise]

Event fires when 'true' if the headstage is connected for the headstage to query has changed

10.104 CSerialPortNet Class Reference

Inheritance diagram for CSerialPortNet:

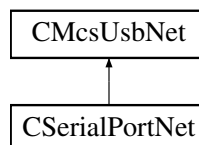

Public Member Functions

- [CSerialPortNet](#) (void)
- void [Send](#) (array< byte > ^ buffer)
- void [Send](#) (String ^ command)
- array< byte > ^ [Receive](#) (void)
- array< byte > ^ [Receive](#) (int length)
- String ^ [ReceiveString](#) (void)
- String ^ [ReceiveString](#) (int length)
- int [GetBytesAvailable](#) (void)

Additional Inherited Members

10.104.1 Constructor & Destructor Documentation

10.104.1.1 CSerialPortNet() `CSerialPortNet (`
`void)`

10.104.2 Member Function Documentation

10.104.2.1 GetBytesAvailable() `int GetBytesAvailable (`
`void)`

10.104.2.2 Receive() [1/2] `array<byte> ^ Receive (`
`int length)`

10.104.2.3 Receive() [2/2] `array<byte> ^ Receive (`
`void)`

10.104.2.4 ReceiveString() [1/2] `String ^ ReceiveString (`
`int length)`

10.104.2.5 ReceiveString() [2/2] `String ^ ReceiveString (`
`void)`


10.104.2.6 Send() [1/2] `void Send (`
`array< byte > ^ buffer)`

10.104.2.7 Send() [2/2] `void Send (`
`String ^ command)`

10.105 CStg200xBasicNet Class Reference

Base class for the Stg200x.

Inheritance diagram for CStg200xBasicNet:

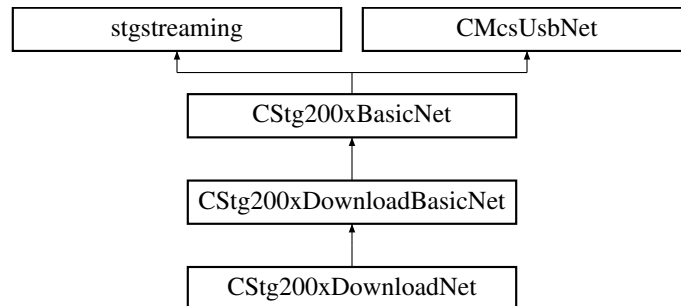

Public Member Functions

- virtual [~CStg200xBasicNet](#) ()
The destructor.
- void [SetOutputRate](#) (uint32_t rate)
Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.
- uint32_t [GetOutputRate](#) ()
Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- void [SendStop](#) (uint32_t triggermap, int options)
Stop some or all triggers of the STG.
- void [GetStgVersionInfo](#) ([Out]String^% SwVersion, [Out]String^% HwVersion)
Queries software and hardware version.
- void [GetAnalogRanges](#) (int channel, [Out]int% URange, [Out]int% IRange)
Gets the range of the analog outputs.
- void [GetAnalogResolution](#) (int channel, [Out]int% URes, [Out]int% IRes)
Gets the resolution of the analog outputs.
- virtual int32_t [GetDACResolution](#) ()
Gets number of bits of the DAC resolution.
- virtual int32_t [GetVoltageRangeInMicroVolt](#) (uint32_t channel)
Gets the Voltage Range of the specified channel in Microvolts.
- virtual int32_t [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)
Gets the Voltage Resolution of the specified channel in Microvolts.
- virtual int32_t [GetCurrentRangeInNanoAmp](#) (uint32_t channel)
Gets the Current Range of the specified channel in Nanoamps.
- virtual int32_t [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)
Gets the Current Resolution of the specified channel in Nanoamps.
- void [GetStgProgramInfo](#) ([Out]bool% IsProgrammed, [Out]System::Runtime::InteropServices::ComTypes::FILETIME% timestamp, [Out]String^% filename, [Out]Guid% guid)
Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.

- void [GetStgProgramInfo](#) ([Out]bool% IsProgrammed, [Out]DateTime% timestamp, [Out]String^% filename, [Out]Guid% guid)
Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.
- void [SetStgProgramInfo](#) (DateTime timestamp, String^ filename, Guid guid)
Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.
- uint32_t [GetAvailableMemory](#) ()
Gets the amount of memory available in the currently selected segment of the STG.
- uint32_t [GetTotalMemory](#) ()
Gets the total amount of memory available on the STG (all segments).
- virtual uint32_t [GetNumberOfAnalogChannels](#) ()
Gets the Number of available analog channels of the device.
- virtual uint32_t [GetNumberOfSyncoutChannels](#) ()
Gets the Number of available syncout channels of the device.
- virtual uint32_t [GetNumberOfTriggerInputs](#) ()
Gets the Number of trigger inputs of the device.
- virtual uint32_t [GetNumberOfHWDACPaths](#) ()
Gets the Number of HW Stimulation DACs of the device.
- virtual uint32_t [GetNumberOfStimulationSourcesPerElectrode](#) ()
Gets the number of stimulation sources (DACs) per electrode.
- virtual void [SetVoltageMode](#) (unsigned int channel)
Sets a channel to voltage mode (STG3008-FA and STG400x only).
- virtual void [SetCurrentMode](#) (unsigned int channel)
Sets a channel to current mode (STG3008-FA and STG400x only).
- virtual void [SetVoltageMode](#) ()
Sets all channels to voltage mode (STG3008-FA and STG400x only).
- virtual void [SetCurrentMode](#) ()
Sets all channels to current mode (STG3008-FA and STG400x only).
- virtual void [SetMeasurementMode](#) (unsigned int channel)
Sets a channel to measurement mode (STG3008-FA).
- virtual void [SetFAAmplification](#) (unsigned int amplification)
- virtual uint32_t [GetFAAmplification](#) ()
- virtual void [SetAutocalibrationDisabled](#) (unsigned int channel, bool disable)
Sets the autocalibration configuration.
- virtual bool [GetAutocalibrationDisabled](#) (unsigned int channel)
Gets the autocalibration configuration.
- virtual void [SetElectrodeMode](#) (uint32_t electrode, array< ElectrodeModeEnumNet >^ mode)
Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t electrode, ElectrodeModeEnumNet mode)
Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode, array< ElectrodeModeEnumNet >^ mode)
Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode, ElectrodeModeEnumNet mode)
Puts an electrode in either automatic or manual mode.
- virtual uint32_t [GetElectrodeMode](#) (uint32_t electrode)
Gets the mode an electrode is in.
- virtual uint32_t [GetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode)
Gets the mode an electrode is in.
- virtual void [SetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodeIndex, array< ElectrodeDacMuxEnumNet >^ dacMux)

Defines the DAC to use for an electrode.

- virtual void [SetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodelIndex, ElectrodeDacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

- virtual void [SetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, ElectrodeDacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

- virtual void [SetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, array< ElectrodeDacMuxEnumNet >^ dacMux)

Defines the DAC to use for an electrode.

- virtual ElectrodeDacMuxEnumNet [GetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodelIndex)

Gets the DAC which is used for an electrode.

- virtual ElectrodeDacMuxEnumNet [GetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex)

Gets the DAC which is used for an electrode.

- virtual void [SetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an electrode.

- virtual void [SetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

- virtual void [SetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

- virtual void [SetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an electrode.

- virtual bool [GetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex)

Gets whether an electrode is enabled or disabled for stimulation.

- virtual bool [GetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex)

Gets whether an electrode is enabled or disabled for stimulation.

- virtual void [SetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an external electrode.

- virtual void [SetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, bool enable)

Enables or disables the stimulation switch for an external electrode.

- virtual bool [GetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex)

Gets whether an electrode is enabled or disabled for stimulation.

- virtual void [SetBlankingEnable](#) (uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

- virtual void [SetBlankingEnable](#) (uint32_t electrode, array< bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

- virtual void [SetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

- virtual void [SetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

- virtual bool [GetBlankingEnable](#) (uint32_t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

- virtual bool [GetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t electrode, bool enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t electrode, array< bool >^ enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode, bool enable)

- Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.*

 - virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)
- Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.*

 - virtual bool [GetEnableAmplifierProtectionSwitch](#) (uint32_t electrode)
- Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.*

 - virtual bool [GetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode)
- Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.*

 - virtual uint32_t [GetNumberOfStimulationElectrodes](#) ()
- template<typename digitalsourceenum >

 - virtual void [SetTriggerSource](#) (unsigned int triggernum, [DigitalSource](#)< digitalsourceenum >^ triggersource, int bitnum_offset)
 - virtual void [SetTriggerSource](#) (unsigned int triggernum, TriggerSourceEnumNet triggersource, int bitnum_offset)
 - virtual void [SetTriggerSource](#) (unsigned int triggernum, TriggerSourceEnumNet triggersource)
 - virtual TriggerSourceEnumNet [GetTriggerSource](#) (unsigned int triggernum)
 - virtual void [SetListmodeIndexRange](#) (unsigned int sideband, unsigned int startIndex, unsigned int endIndex, unsigned int mode)
 - virtual void [GetListmodeIndexRange](#) (unsigned int sideband, unsigned int &startIndex, unsigned int &endIndex, unsigned int &mode)
 - virtual void [SetListmodeTriggerSource](#) (unsigned int sideband, TriggerSourceEnumNet triggersource, int bitnumOffset)
 - virtual void [SetListmodeTriggerSource](#) (unsigned int sideband, TriggerSourceEnumNet triggersource)
 - virtual TriggerSourceEnumNet [GetListmodeTriggerSource](#) (unsigned int sideband)
 - virtual void [ListModeSendStart](#) (unsigned int sidebandMask)
 - virtual void [ListModeSendStop](#) (unsigned int sidebandMask)
 - virtual void [SetHeadstage](#) (unsigned int headstage)
 - virtual uint32_t [GetHeadstage](#) ()
 - virtual void [SetDacAmplificationFactor](#) (uint32_t DacNumber, double Factor)
- Set the amplification factor for a DAC.*

 - virtual double [GetDacAmplificationFactor](#) (uint32_t DacNumber)
- Get the amplification factor for a DAC.*

 - virtual void [SetDigoutMode](#) (Stg200xDigoutModeEnumNet digoutMode)
- Sets the operation mode of the digital output port, can be Monitor, Manual or SyncOut*

 - virtual Stg200xDigoutModeEnumNet [GetDigoutMode](#) ()
- Gets the operation mode of the digital output port, can be Monitor, Manual or SyncOut*

 - virtual void [SetDigoutValue](#) (uint32_t digoutValue)
- Sets the Value on the digital output port when in manual mode.*

 - virtual uint32_t [GetDigoutValue](#) ()
- Gets the Value on the digital output port.*

 - virtual uint32_t [GetDiginValue](#) ()
- Gets the Value on the digital input port.*

 - virtual void [SetSyncoutMap](#) (uint32_t channel, uint32_t syncoutMap)
- Sets the mapping between external syncout outputs and internal syncout channels.*

 - virtual uint32_t [GetSyncoutMap](#) (uint32_t channel)
- Gets the mapping between external syncout outputs and internal syncout channels.*

Additional Inherited Members

10.105.1 Detailed Description

Base class for the Stg200x.

From this class all STG related classes are derived: UsbNetDll::CStg200xDownloadBasicNet UsbNetDll::CStg200xDownloadNet for [Download Mode](#) and UsbNetDll::CStg200xStreamingNet for [Streaming Mode](#).

[CStg200xBasicNet](#) is the base class to control MCS STG device.

10.105.2 Constructor & Destructor Documentation

10.105.2.1 `~CStg200xBasicNet()` `virtual ~CStg200xBasicNet () [virtual]`

The destructor.

10.105.3 Member Function Documentation

10.105.3.1 `GetAnalogRanges()` `void GetAnalogRanges (`
`int channel,`
`[Out] int% URange,`
`[Out] int% IRange)`

Gets the range of the analog outputs.

Parameters

<i>channel</i>	The channel which is queried.
<i>URange</i>	The Voltage range in mV.
<i>IRange</i>	The Current range in uA.

10.105.3.2 `GetAnalogResolution()` `void GetAnalogResolution (`
`int channel,`

```
[Out] int% URes,  
[Out] int% IRes )
```

Gets the resolution of the analog outputs.

Parameters

<i>channel</i>	The channel which is queried.
<i>URes</i>	The Voltage resolution in mV.
<i>IRes</i>	The Current resolution in uA.

10.105.3.3 GetAutocalibrationDisabled() `virtual bool GetAutocalibrationDisabled (unsigned int channel) [virtual]`

Gets the autocalibration configuration.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

`true` if autocalibration is disabled.

10.105.3.4 GetAvailableMemory() `uint32_t GetAvailableMemory ()`

Gets the amount of memory available in the currently selected segment of the STG.

Returns

The memory available in the currently selected segment in bytes.

10.105.3.5 GetBlankingEnable() `[1/2] virtual bool GetBlankingEnable (uint32_t electrode) [virtual]`

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if blanking is enabled while stimulation is in progress.

10.105.3.6 GetBlankingEnable() [2/2] `virtual bool GetBlankingEnable (`
`uint32_t Scu_HS,`
`uint32_t electrode) [virtual]`

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if blanking is enabled while stimulation is in progress.

10.105.3.7 GetCurrentRangeInNanoAmp() `virtual int32_t GetCurrentRangeInNanoAmp (`
`uint32_t channel) [virtual]`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

10.105.3.8 GetCurrentResolutionInNanoAmp() `virtual int32_t GetCurrentResolutionInNanoAmp (`
`uint32_t channel) [virtual]`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

10.105.3.9 GetDacAmplificationFactor() `virtual double GetDacAmplificationFactor (uint32_t DacNumber) [virtual]`

Get the amplification factor for a DAC.

Parameters

<i>DacNumber</i>	The number of the DAC.
------------------	------------------------

Returns

the amplification factor for the DAC queried, range is from -1.99999 to +1.99999.

10.105.3.10 GetDACResolution() `virtual int32_t GetDACResolution () [virtual]`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

10.105.3.11 GetDiginValue() `virtual uint32_t GetDiginValue () [virtual]`

Gets the Value on the digital input port.

Returns

The current value on the digital inputs.

10.105.3.12 GetDigoutMode() `virtual Stg200xDigoutModeEnumNet GetDigoutMode () [virtual]`

Gets the operation mode of the digital output port, can be Monitor, Manual or SyncOut

Returns

The current operation mode.

10.105.3.13 GetDigoutValue() `virtual uint32_t GetDigoutValue () [virtual]`

Gets the Value on the digital output port.

Returns

The current value on the digital outputs.

10.105.3.14 GetElectrodeDacMux() [1/2] `virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex) [virtual]`

Gets the DAC which is used for an electrode.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

10.105.3.15 GetElectrodeDacMux() [2/2] `virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex) [virtual]`

Gets the DAC which is used for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

10.105.3.16 GetElectrodeEnable() [1/2] virtual bool GetElectrodeEnable (
 uint32_t *electrode*,
 uint32_t *listmodeIndex*) [virtual]

Gets whether an electrode is enabled or disabled for stimulation.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

10.105.3.17 GetElectrodeEnable() [2/2] virtual bool GetElectrodeEnable (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 uint32_t *listmodeIndex*) [virtual]

Gets whether an electrode is enabled or disabled for stimulation.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

10.105.3.18 GetElectrodeMode() [1/2] `virtual uint32_t GetElectrodeMode (uint32_t electrode) [virtual]`

Gets the mode an electrode is in.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

10.105.3.19 GetElectrodeMode() [2/2] `virtual uint32_t GetElectrodeMode (uint32_t Scu_HS, uint32_t electrode) [virtual]`

Gets the mode an electrode is in.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

10.105.3.20 GetEnableAmplifierProtectionSwitch() [1/2] `virtual bool GetEnableAmplifierProtectionSwitch (`
`uint32_t electrode) [virtual]`

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

10.105.3.21 GetEnableAmplifierProtectionSwitch() [2/2] `virtual bool GetEnableAmplifierProtectionSwitch (`
`uint32_t Scu_HS,`
`uint32_t electrode) [virtual]`

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

10.105.3.22 GetExternalElectrodeEnable() `virtual bool GetExternalElectrodeEnable (`
`uint32_t electrode,`
`uint32_t listmodeIndex) [virtual]`

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

10.105.3.23 GetFAAmplification() `virtual uint32_t GetFAAmplification () [virtual]`

10.105.3.24 GetHeadstage() `virtual uint32_t GetHeadstage () [virtual]`

10.105.3.25 GetListmodeIndexRange() `virtual void GetListmodeIndexRange (`
 `unsigned int sideband,`
 `unsigned int & startIndex,`
 `unsigned int & endIndex,`
 `unsigned int & mode) [virtual]`

10.105.3.26 GetListmodeTriggerSource() `virtual TriggerSourceEnumNet GetListmodeTriggerSource (`
 `unsigned int sideband) [virtual]`

10.105.3.27 GetNumberOfAnalogChannels() `virtual uint32_t GetNumberOfAnalogChannels () [virtual]`

Gets the Number of available analog channels of the device.

Returns

The number of analog channels.

10.105.3.28 GetNumberOfHWDACPaths() `virtual uint32_t GetNumberOfHWDACPaths () [virtual]`

Gets the Number of HW Stimulation DACs of the device.

Returns

The number of independent HW Stimulation outputs.

10.105.3.29 GetNumberOfStimulationElectrodes() virtual uint32_t GetNumberOfStimulationElectrodes () [virtual]

10.105.3.30 GetNumberOfStimulationSourcesPerElectrode() virtual uint32_t GetNumberOfStimulationSourcesPerElectrode () [virtual]

Gets the number of stimulation sources (DACs) per electrode.

Returns

The number of stimulation sources (DACs) per electrode.

10.105.3.31 GetNumberOfSyncoutChannels() virtual uint32_t GetNumberOfSyncoutChannels () [virtual]

Gets the Number of available syncout channels of the device.

Returns

The number of analog channels.

10.105.3.32 GetNumberOfTriggerInputs() virtual uint32_t GetNumberOfTriggerInputs () [virtual]

Gets the Number of trigger inputs of the device.

Returns

The number of trigger inputs.

10.105.3.33 GetOutputRate() uint32_t GetOutputRate ()

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Returns

Returns the current output rate in Hz.

10.105.3.34 GetStgProgramInfo() [1/2] void GetStgProgramInfo ([Out] bool% *IsProgrammed*, [Out] DateTime% *timestamp*, [Out] String^% *filename*, [Out] Guid% *guid*)

Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.

Parameters

<i>IsProgrammed</i>	Flag wether download information is valid.
<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.
<i>guid</i>	A GUID.

10.105.3.35 GetStgProgramInfo() [2/2] `void GetStgProgramInfo (`
`[Out] bool% IsProgrammed,`
`[Out] System::Runtime::InteropServices::ComTypes::FILETIME% timestamp,`
`[Out] String^% filename,`
`[Out] Guid% guid)`

Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.

Parameters

<i>IsProgrammed</i>	Flag wether download information is valid.
<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.

10.105.3.36 GetStgVersionInfo() `void GetStgVersionInfo (`
`[Out] String^% SwVersion,`
`[Out] String^% HwVersion)`

Queries software and hardware version.

Parameters

<i>SwVersion</i>	The current Software Version of the STG.
<i>HwVersion</i>	The Hardware Revision of the STG.

10.105.3.37 GetSyncoutMap() `virtual uint32_t GetSyncoutMap (`
`uint32_t channel) [virtual]`

Gets the mapping between external syncout outputs and internal syncout channels.

Parameters

<i>channel</i>	The external syncout output channel number (zero based).
----------------	--

Returns

The bitmap of internal syncout channels mapped to channel.

10.105.3.38 GetTotalMemory() `uint32_t GetTotalMemory ()`

Gets the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

10.105.3.39 GetTriggerSource() `virtual TriggerSourceEnumNet GetTriggerSource (unsigned int triggernum) [virtual]`**10.105.3.40 GetVoltageRangeInMicroVolt()** `virtual int32_t GetVoltageRangeInMicroVolt (uint32_t channel) [virtual]`

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

10.105.3.41 GetVoltageResolutionInMicroVolt() `virtual int32_t GetVoltageResolutionInMicroVolt (uint32_t channel) [virtual]`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

10.105.3.42 ListModeSendStart() virtual void ListModeSendStart (
 unsigned int *sidebandMask*) [virtual]

10.105.3.43 ListModeSendStop() virtual void ListModeSendStop (
 unsigned int *sidebandMask*) [virtual]

10.105.3.44 SendStart() void SendStart (
 uint32_t *triggermap*)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

10.105.3.45 SendStop() [1/2] void SendStop (
 uint32_t *triggermap*)

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

10.105.3.46 SendStop() [2/2] void SendStop (
 uint32_t *triggermap*,
 int *options*)

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
<i>options</i>	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout associated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active..

10.105.3.47 SetAutocalibrationDisabled() `virtual void SetAutocalibrationDisabled (unsigned int channel, bool disable) [virtual]`

Sets the autocalibration configuration.

Parameters

<i>channel</i>	The channel number.
<i>disable</i>	true if autocalibration is to be disabled.

10.105.3.48 SetBlankingEnable() [1/4] `virtual void SetBlankingEnable (uint32_t electrode, array< bool >^ enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.49 SetBlankingEnable() [2/4] `virtual void SetBlankingEnable (uint32_t electrode, bool enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.50 SetBlankingEnable() [3/4] `virtual void SetBlankingEnable (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.51 SetBlankingEnable() [4/4] `virtual void SetBlankingEnable (`
 `uint32_t Scu_HS,`
 `uint32_t electrode,`
 `bool enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.52 SetCurrentMode() [1/2] `virtual void SetCurrentMode () [virtual]`

Sets all channels to current mode (STG3008-FA and STG400x only).

10.105.3.53 SetCurrentMode() [2/2] `virtual void SetCurrentMode (`
 `unsigned int channel) [virtual]`

Sets a channel to current mode (STG3008-FA and STG400x only).

Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

10.105.3.54 SetDacAmplificationFactor() virtual void SetDacAmplificationFactor (
 uint32_t *DacNumber*,
 double *Factor*) [virtual]

Set the amplification factor for a DAC.

Parameters

<i>DacNumber</i>	The number of the DAC.
<i>Factor</i>	the amplification factor for that DAC, range is from -1.99999 to +1.99999.

10.105.3.55 SetDigoutMode() virtual void SetDigoutMode (
 Stg200xDigoutModeEnumNet *digoutMode*) [virtual]

Sets the operation mode of the digital output port, can be Monitor, Manual or SyncOut

Parameters

<i>digoutMode</i>	The new operation mode.
-------------------	-------------------------

10.105.3.56 SetDigoutValue() virtual void SetDigoutValue (
 uint32_t *digoutValue*) [virtual]

Sets the Value on the digital output port when in manual mode.

Parameters

<i>digoutValue</i>	The new value on the digital outputs.
--------------------	---------------------------------------

10.105.3.57 SetElectrodeDacMux() [1/4] virtual void SetElectrodeDacMux (
 uint32_t *electrode*,
 uint32_t *listmodeIndex*,
 array< ElectrodeDacMuxEnumNet >^ *dacMux*) [virtual]

Defines the DAC to use for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

10.105.3.58 SetElectrodeDacMux() [2/4] virtual void SetElectrodeDacMux (
 uint32_t *electrode*,
 uint32_t *listmodeIndex*,
 ElectrodeDacMuxEnumNet *dacMux*) [virtual]

Defines the DAC to use for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

10.105.3.59 SetElectrodeDacMux() [3/4] virtual void SetElectrodeDacMux (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 uint32_t *listmodeIndex*,
 array< ElectrodeDacMuxEnumNet >^ *dacMux*) [virtual]

Defines the DAC to use for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

10.105.3.60 SetElectrodeDacMux() [4/4] `virtual void SetElectrodeDacMux (`
`uint32_t Scu_HS,`
`uint32_t electrode,`
`uint32_t listmodeIndex,`
`ElectrodeDacMuxEnumNet dacMux) [virtual]`

Defines the DAC to use for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

10.105.3.61 SetElectrodeEnable() [1/4] `virtual void SetElectrodeEnable (`
`uint32_t electrode,`


```
uint32_t listmodeIndex,
array< bool >^ enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.62 SetElectrodeEnable() [2/4] virtual void SetElectrodeEnable (

```
uint32_t electrode,
uint32_t listmodeIndex,
bool enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.63 SetElectrodeEnable() [3/4] virtual void SetElectrodeEnable (

```
uint32_t Scu_HS,
uint32_t electrode,
uint32_t listmodeIndex,
array< bool >^ enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.64 SetElectrodeEnable() [4/4] `virtual void SetElectrodeEnable (`
 `uint32_t Scu_HS,`
 `uint32_t electrode,`
 `uint32_t listmodeIndex,`
 `bool enable) [virtual]`

Enables or disables the stimulation switch for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.65 SetElectrodeMode() [1/4] virtual void SetElectrodeMode (
 uint32_t *electrode*,
 array< ElectrodeModeEnumNet >^ *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

10.105.3.66 SetElectrodeMode() [2/4] virtual void SetElectrodeMode (
 uint32_t *electrode*,
 ElectrodeModeEnumNet *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>mode</i>	0 for automatic and 3 for manual mode.
-------------	--

10.105.3.67 SetElectrodeMode() [3/4] virtual void SetElectrodeMode (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 array< ElectrodeModeEnumNet >^ *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

10.105.3.68 SetElectrodeMode() [4/4] virtual void SetElectrodeMode (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 ElectrodeModeEnumNet *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>mode</i>	0 for automatic and 3 for manual mode.
-------------	--

10.105.3.69 SetEnableAmplifierProtectionSwitch() [1/4] virtual void SetEnableAmplifierProtection↔
 Switch (

```
uint32_t electrode,
array< bool >^ enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.70 SetEnableAmplifierProtectionSwitch() [2/4] virtual void SetEnableAmplifierProtection↔
Switch (

```
uint32_t electrode,
bool enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.71 SetEnableAmplifierProtectionSwitch() [3/4] virtual void SetEnableAmplifierProtection↔
Switch (

```
uint32_t Scu_HS,
uint32_t electrode,
array< bool >^ enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.72 SetEnableAmplifierProtectionSwitch() [4/4] virtual void SetEnableAmplifierProtection↔
Switch (

```
uint32_t Scu_HS,
uint32_t electrode,
bool enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

10.105.3.73 SetExternalElectrodeEnable() [1/2] virtual void SetExternalElectrodeEnable (
uint32_t electrode,
uint32_t listmodeIndex,
array< bool >^ enable) [virtual]

Enables or disables the stimulation switch for an external electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.74 SetExternalElectrodeEnable() [2/2] virtual void SetExternalElectrodeEnable (
uint32_t electrode,
uint32_t listmodeIndex,
bool enable) [virtual]

Enables or disables the stimulation switch for an external electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

10.105.3.75 SetFAAmplification() virtual void SetFAAmplification (
 unsigned int *amplification*) [virtual]

10.105.3.76 SetHeadstage() virtual void SetHeadstage (
 unsigned int *headstage*) [virtual]

10.105.3.77 SetListmodeIndexRange() virtual void SetListmodeIndexRange (
 unsigned int *sideband*,
 unsigned int *startIndex*,
 unsigned int *endIndex*,
 unsigned int *mode*) [virtual]

10.105.3.78 SetListmodeTriggerSource() [1/2] virtual void SetListmodeTriggerSource (
 unsigned int *sideband*,
 TriggerSourceEnumNet *triggersource*) [virtual]

10.105.3.79 SetListmodeTriggerSource() [2/2] virtual void SetListmodeTriggerSource (
 unsigned int *sideband*,
 TriggerSourceEnumNet *triggersource*,
 int *bitnumOffset*) [virtual]

10.105.3.80 SetMeasurementMode() virtual void SetMeasurementMode (
 unsigned int *channel*) [virtual]

Sets a channel to measurement mode (STG3008-FA).

Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

10.105.3.81 SetOutputRate() `void SetOutputRate (`
`uint32_t rate)`

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Parameters

<i>rate</i>	The new output rate in Hz.
-------------	----------------------------

10.105.3.82 SetStgProgramInfo() `void SetStgProgramInfo (`
`DateTime timestamp,`
`String^ filename,`
`Guid guid)`

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

Parameters

<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.

10.105.3.83 SetSyncoutMap() `virtual void SetSyncoutMap (`
`uint32_t channel,`
`uint32_t syncoutMap) [virtual]`

Sets the mapping between external syncout outputs and internal syncout channels.

Parameters

<i>channel</i>	The external syncout output channel number (zero based).
<i>syncoutMap</i>	A bitmap of internal syncout channels to map to channel.

10.105.3.84 SetTriggerSource() [1/3] `virtual void SetTriggerSource (`
`unsigned int triggernum,`


```
DigitalSource< digitalsourceenum >^ triggersource,
int bitnum_offset ) [virtual]
```

10.105.3.85 SetTriggerSource() [2/3] virtual void SetTriggerSource (unsigned int *triggernum*, TriggerSourceEnumNet *triggersource*) [virtual]

10.105.3.86 SetTriggerSource() [3/3] virtual void SetTriggerSource (unsigned int *triggernum*, TriggerSourceEnumNet *triggersource*, int *bitnum_offset*) [virtual]

10.105.3.87 SetVoltageMode() [1/2] virtual void SetVoltageMode () [virtual]

Sets all channels to voltage mode (STG3008-FA and STG400x only).

10.105.3.88 SetVoltageMode() [2/2] virtual void SetVoltageMode (unsigned int *channel*) [virtual]

Sets a channel to voltage mode (STG3008-FA and STG400x only).


Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

10.106 CStg200xDownloadBasicNet Class Reference

[CStg200xDownloadBasicNet](#) is the base class to control the download mode of the MCS STG device.

Inheritance diagram for CStg200xDownloadBasicNet:

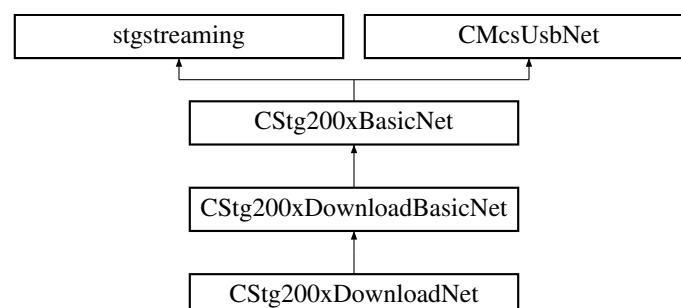

Public Member Functions

- virtual void [SetupTrigger](#) (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)
Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- virtual void [SetupTriggerSingle](#) (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)
Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void [GetTrigger](#) ([Out] array< uint32_t >^% channelmap, [Out] array< uint32_t >^% syncoutmap, [Out] array< uint32_t >^% repeat)
Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void [GetSweepCount](#) ([Out] array< uint32_t >^% sweeps, [Out] array< uint32_t >^% triggers)
Get the sweep and trigger count of the STG.
 - The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
 - The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in [CStg200xDownloadBasicNet::SetupTrigger](#).
- void [ForceStatusEvent](#) ()
Force a status event.
- void [ResetStatus](#) (uint32_t triggermap)
Reset the status flag.
- uint32_t [GetMemoryUsageDAC](#) (uint32_t Channel)
Queries the memory usage of the current segment and selected analog DAC channel.
- uint32_t [GetMemoryUsageSyncout](#) (uint32_t Channel)
Queries the memory usage of the current segment and selected syncout channel.
- virtual void [ClearSyncData](#) (uint32_t channel)
Delete a SyncOut pattern for a channel from STG memory.
- virtual void [SendSyncData](#) (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)
Uploads sync output data to the STG.
Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.
Each datapoint is represented by an integer value and can be either 0 or 1.
The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the [STG200x_SendSyncData32\(\)](#) call instead.
- virtual void [ClearChannelData](#) (uint32_t channel)
Delete a stimulus pattern for a channel from STG memory
- virtual void [SendChannelData](#) (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)
Uploads analog data (stimulus patterns) to the STG.
Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.
Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.
The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.
- virtual void [EnableAutoReset](#) ()
Enable AutoReset of the STG Status.
- virtual void [DisableAutoReset](#) ()
Disable AutoReset of the STG Status.
- virtual void [SetupRetriggerMode](#) (int8_t trigger, RetriggerActionEnumNet same_trigger, RetriggerActionEnumNet other_trigger)
Define the action on triggers while the STG is running.
The STG has three options how to handle a successive trigger while a trigger is active.
 - stop this trigger (default action)
 - restart this trigger
 - ignore the signal
- virtual void [SetupRetriggerMode](#) (RetriggerActionEnumNet same_trigger, RetriggerActionEnumNet other_trigger)

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- *stop this trigger (default action)*
- *restart this trigger*
- *ignore the signal*

Properties

- [CStimulusFunctionNet](#)[^] [Stimulus](#) [get]

Additional Inherited Members

10.106.1 Detailed Description

[CStg200xDownloadBasicNet](#) is the base class to control the download mode of the MCS STG device.

10.106.2 Member Function Documentation

10.106.2.1 ClearChannelData() `virtual void ClearChannelData (
 uint32_t channel) [virtual]`

Delete a stimulus pattern for a channel from STG memory

Parameters

<i>channel</i>	Specifies the channel to clear.
----------------	---------------------------------

10.106.2.2 ClearSyncData() `virtual void ClearSyncData (
 uint32_t channel) [virtual]`

Delete a SyncOut pattern for a channel from STG memory.

Parameters

<i>channel</i>	Specifies the syncout channel to clear.
----------------	---

10.106.2.3 DisableAutoReset() `virtual void DisableAutoReset () [virtual]`

Disable AutoReset of the STG Status.

If autoreset is disabled, the STG status switches to FINISHED after the defined number of sweeps is finished. To switch back to the IDLE status, use CStg200xDownload::ResetStatus()

10.106.2.4 EnableAutoReset() `virtual void EnableAutoReset () [virtual]`

Enable AutoReset of the STG Status.

This is the default on power up. If autoreset is enabled, the STG status switches to FINISHED only for one poll cycle after this, it switches to IDLE automatically.

10.106.2.5 ForceStatusEvent() `void ForceStatusEvent ()`

Force a status event.

Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

10.106.2.6 GetMemoryUsageDAC() `uint32_t GetMemoryUsageDAC (uint32_t Channel)`

Queries the memory usage of the current segment and selected analog DAC channel.

The currently used memory is reported for the requested channel.

Parameters

<i>Channel</i>	channel for the amount of interested usage.
----------------	---

Returns

Returns the usage in STG memory.

10.106.2.7 GetMemoryUsageSyncout() `uint32_t GetMemoryUsageSyncout (uint32_t Channel)`

Queries the memory usage of the current segment and selected syncout channel.

The currently used memory is reported for the requested channel.

Parameters

<i>Channel</i>	channel for the amount of interested usage.
----------------	---

Returns

Returns the usage in STG memory.

10.106.2.8 GetSweepCount() `void GetSweepCount (`
 `[Out] array< uint32_t > ^% sweeps,`
 `[Out] array< uint32_t > ^% triggers)`

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in [CStg200xDownloadBasicNet::SetupTrigger](#).

Parameters

<i>sweeps</i>	on return contains the number of sweeps for each trigger.
<i>triggers</i>	on return contains the number of trigger events seen for each trigger.

10.106.2.9 GetTrigger() `void GetTrigger (`
 `[Out] array< uint32_t > ^% channelmap,`
 `[Out] array< uint32_t > ^% syncoutmap,`
 `[Out] array< uint32_t > ^% repeat)`

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

10.106.2.10 ResetStatus() `void ResetStatus (`
 `uint32_t triggermap)`

Reset the status flag.

Parameters

<i>triggermap</i>	bitmap of trigger for which to reset the status.
-------------------	--

10.106.2.11 SendChannelData() virtual void SendChannelData (
 uint32_t *channel*,
 array< uint16_t >^ *pData*,
 array< uint64_t >^ *tData*) [virtual]

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Parameters

<i>channel</i>	Specifies the channel to append the data to.
<i>pData</i>	A list of datapoints.
<i>tData</i>	A list of durations as int64_t. The time is given in units of μ s.

10.106.2.12 SendSyncData() virtual void SendSyncData (
 uint32_t *channel*,
 array< uint16_t >^ *pData*,
 array< uint64_t >^ *tData*) [virtual]

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the STG200x_SendSyncData32() call instead.

Parameters

<i>channel</i>	Specifies the sync output channel to append the data to.
<i>pData</i>	A list of datapoints.
<i>tData</i>	A list of durations as int64_t. The time is given in units of μ s.

10.106.2.13 SetupRetriggerMode() [1/2] `virtual void SetupRetriggerMode (`
`int8_t trigger,`
`RetriggerActionEnumNet same_trigger,`
`RetriggerActionEnumNet other_trigger) [virtual]`

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Parameters

<i>trigger</i>	The trigger to change.
<i>same_trigger</i>	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
<i>other_trigger</i>	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

10.106.2.14 SetupRetriggerMode() [2/2] `virtual void SetupRetriggerMode (`
`RetriggerActionEnumNet same_trigger,`
`RetriggerActionEnumNet other_trigger) [virtual]`

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Parameters

<i>same_trigger</i>	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
<i>other_trigger</i>	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

10.106.2.15 SetupTrigger() `virtual void SetupTrigger (`


```
uint32_t first_trigger,
array< uint32_t >^ channelmap,
array< uint32_t >^ syncoutmap,
array< uint32_t >^ repeat ) [virtual]
```

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>first_trigger</i>	The number of the first trigger to change.
----------------------	--

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

10.106.2.16 SetupTriggerSingle() virtual void SetupTriggerSingle (

```
uint32_t trigger,
uint32_t channelmap,
uint32_t syncoutmap,
uint32_t repeat ) [virtual]
```

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>trigger</i>	The trigger to change.
----------------	------------------------

Parameters

<i>channelmap</i>	A bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	A bitmap of syncouts that belong to this trigger.
<i>repeat</i>	The number of times this trigger should be repeated.


10.106.3 Property Documentation

10.106.3.1 Stimulus [CStimulusFunctionNet](#)[^] Stimulus [get]

10.107 CStg200xDownloadNet Class Reference

Main class for the STG download mode This class implements the STG download mode interface.

Inheritance diagram for CStg200xDownloadNet:

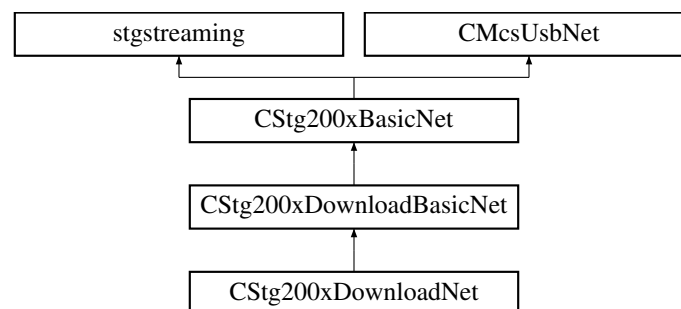

Public Member Functions

- [CStg200xDownloadNet](#) ()
Use this constructor if you do not want to use the status callback.
- [CStg200xDownloadNet](#) ([OnStgPollStatus](#)[^] pollStatus)
Use this constructor if you want to use the status callback.
- [~CStg200xDownloadNet](#) ()
- void [PrepareAndSendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType)
Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.
- void [PrepareAndAppendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType)
Prepare and append data to a given channel on the STG.
- void [ClearChannel_PrepareAndSendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType, bool doClear)
Prepare and append data to a given channel on the STG.
- void [SendSegmentDefine](#) (array< uint32_t >[^] segment_list)

Defines the segment memory layout of the STG.

- void [SendSegmentStart](#) (uint32_t triggermap, uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)

Switchs segment and starts trigger.

- void [SendSegmentSelect](#) (uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)

Switchs segment.

- void [EnableMultiFileMode](#) (uint32_t submode)

Enable the Multi-File mode of the STG.

- void [DisableMultiFileMode](#) ()

Disable the Multi-File mode of the STG

- [StgStatusNet](#) ^ [QueryTriggerstatus](#) ()
- void [SetOutputMap](#) (uint32_t ChannelLayout[])
- int32_t [GetModuleTemp](#) (unsigned int channel)
- uint32_t [GetModuleCurrent](#) (unsigned int channel)

Events

- [OnStgPollStatus](#) ^ [Stg200xPollStatusEvent](#) [add, remove, raise]
- [OnMwPollStatus](#) ^ [MwPollStatusEvent](#) [add, remove, raise]

Additional Inherited Members

10.107.1 Detailed Description

Main class for the STG download mode This class implements the STG download mode interface.

10.107.2 Constructor & Destructor Documentation

10.107.2.1 CStg200xDownloadNet() [1/2] [CStg200xDownloadNet](#) ()

Use this constructor if you do not want to use the status callback.

10.107.2.2 CStg200xDownloadNet() [2/2] [CStg200xDownloadNet](#) ([OnStgPollStatus](#) ^ *pollStatus*)

Use this constructor if you want to use the status callback.

10.107.2.3 ~CStg200xDownloadNet() [~CStg200xDownloadNet](#) ()

10.107.3 Member Function Documentation

10.107.3.1 ClearChannel_PrepareAndSendData() `void ClearChannel_PrepareAndSendData (`
 `uint32_t channel,`
 `array< int32_t >^ amplitude,`
 `array< uint64_t >^ duration,`
 `STG_DestinationEnumNet destType,`
 `bool doClear)`

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

10.107.3.2 DisableMultiFileMode() `void DisableMultiFileMode ()`

Disable the Multi-File mode of the STG

Switch the STG back to normal mode. In this mode, trigger inputs are assigned to channels, not to segments.

10.107.3.3 EnableMultiFileMode() `void EnableMultiFileMode (`
`uint32_t submode)`

Enable the Multi-File mode of the STG.

In Multi-File mode, the trigger inputs switch between segments. To use this mode, define up to as many segments as trigger inputs are available and fill each segment with a stimulus pattern.

Now a trigger on trigger input 1 switches the STG to the first segment and starts all triggers in this segment. Likewise, a trigger on trigger input 2, 3 and 4 selects the respective segment and start all triggers in this segment. So the Multi-File Mode can be used to predefine up to four different stimuli which can be selected without the need for a computer connection.

Parameters

<i>submode</i>	The submode. Submode 0 is regular Multi-File mode as described above, submode 1 is extended Multi-File mode, where the segment is selected based on the digital pattern on the digital inputs. In this mode, 256 different segments can be defined and used.
----------------	--

10.107.3.4 GetModuleCurrent() `uint32_t GetModuleCurrent (`
`unsigned int channel)`

10.107.3.5 GetModuleTemp() `int32_t GetModuleTemp (`
`unsigned int channel)`

10.107.3.6 PrepareAndAppendData() `void PrepareAndAppendData (`
`uint32_t channel,`
`array< int32_t >^ amplitude,`
`array< uint64_t >^ duration,`
`STG_DestinationEnumNet destType)`

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
------------------	---

Parameters

<i>duration</i>	A list of durations in units of μs .
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

10.107.3.7 PrepareAndSendData() `void PrepareAndSendData (`
`uint32_t channel,`
`array< int32_t >^ amplitude,`
`array< uint64_t >^ duration,`
`STG_DestinationEnumNet destType)`

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μV , thus the possible range is $\pm 2000 \text{ V}$. When using current stimulation, the values are in multiple of 1 nA , this the possible range is $\pm 2000 \text{ mA}$.

The duration is given as a list of 64 bit integers. Durations are given in units of μs . The STG has a resolution of 20 μs .

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such a block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
------------------	---

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

10.107.3.8 QueryTriggerstatus() `StgStatusNet ^ QueryTriggerstatus ()`**10.107.3.9 SendSegmentDefine()** `void SendSegmentDefine (array< uint32_t >^ segment_list)`

Defines the segment memory layout of the STG.

On reset, the STG has one segment containing all available memory.

With this command, the STG memory can be devided into several segments. Each segment can be filled with stimulus data.

Parameters

<i>segment_list</i>	The List of memory sizes (one per segment).
---------------------	---

10.107.3.10 SendSegmentSelect() `void SendSegmentSelect (uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)`

Switchs segment.

Parameters

<i>segment</i>	The number of the segment to select.
----------------	--------------------------------------

Parameters

<i>segmentflags</i>	A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.
---------------------	---

10.107.3.11 SendSegmentStart() `void SendSegmentStart (`
 `uint32_t triggermap,`
 `uint32_t segment,`
 `Stg200xSegmentFlagsEnumNet segmentflags)`

Switchs segment and starts trigger.

Parameters

<i>triggermap</i>	A bitmap of triggers that will be started.
-------------------	--

Parameters

<i>segment</i>	The number of the segment to select.
----------------	--------------------------------------

Parameters

<i>segmentflags</i>	A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.
---------------------	---

10.107.3.12 SetOutputMap() `void SetOutputMap (`
 `uint32_t ChannelLayout [])`


10.107.4 Event Documentation

10.107.4.1 MwPollStatusEvent [OnMwPollStatus](#)[^] MwPollStatusEvent [add], [remove], [raise]

10.107.4.2 Stg200xPollStatusEvent [OnStgPollStatus](#)[^] Stg200xPollStatusEvent [add], [remove], [raise]

10.108 CStimulusFunctionNet Class Reference

Inheritance diagram for CStimulusFunctionNet:

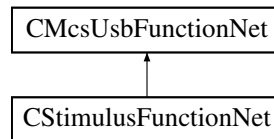

Classes

- class [SidebandData](#)
- class [StimulusDeviceDataAndUnrolledData](#)

Public Member Functions

- [CStimulusFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] stimulusFunctionPointerContainer)
- [CStimulusFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [StartPoll](#) ()
Starts the interrupt fetching thread and delivers events
- void [StopPoll](#) ()
Stops the interrupt fetching thread and delivers events
- void [ForceStatusEvent](#) ()
Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- void [SendStop](#) (uint32_t triggermap, int options)
Stop some or all triggers of the STG.
- void [ClearChannelData](#) (int channel)
Delete a Stimulus Pattern from STG memory
- void [ClearSyncData](#) (int channel)
Delete a Syncout Pattern from STG memory
- void [PrepareAndSendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType)
Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.
- void [PrepareAndAppendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType)
Prepare and append data to a given channel on the STG.
- void [ClearChannel_PrepareAndSendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType, bool doClear)
- [StimulusDeviceDataAndUnrolledData](#)[^] [PrepareData](#) (int channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType)
- void [SendPreparedData](#) (int channel, [StimulusDeviceDataAndUnrolledData](#)[^] device_data_and_unrolled, STG_DestinationEnumNet destType)
- [SidebandData](#)[^] [CreateSideband](#) (array< int32_t >[^] StimulusActive, array< int32_t >[^] Syncout, array< uint64_t >[^] Duration, uint32_t Bit0Time, uint32_t Bit3Time, uint32_t Bit4Time)

- Creates the Sideband Channel for the MEA2100 device.*

 - void [ClearMultiplexedData](#) ()
- Clears the Stimulation Memory in the STG device.*

 - void [SendMultiplexedData](#) (array< uint16_t >^ data)

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.
- int [GetMultiplexedDataChannelsInBlock](#) ()

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData
- int [GetDACResolution](#) ()

Gets number of bits of the DAC resolution.
- int [GetVoltageRangeInMicroVolt](#) (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.
- int [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.
- int [GetCurrentRangeInNanoAmp](#) (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.
- int [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.
- void [SetupTrigger](#) (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void [SetupTriggerSingle](#) (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- uint32_t [GetTotalMemory](#) ()

Get the total amount of memory available on the STG (all segments).
- uint32_t [GetAvailableMemory](#) ()

Get the amount of memory available in the currently selected segment of the STG.
- int [GetNumberOfAnalogChannels](#) ()

Get the number of STG channels.

Events

- [OnStgPollStatus](#)[^] [PollStatusEvent](#)

Additional Inherited Members

10.108.1 Constructor & Destructor Documentation

10.108.1.1 CStimulusFunctionNet() [1/2] [CStimulusFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] stimulusFunctionPointerContainer)

10.108.1.2 CStimulusFunctionNet() [2/2] [CStimulusFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

10.108.2 Member Function Documentation

10.108.2.1 ClearChannel_PrepareAndSendData() `void ClearChannel_PrepareAndSendData (`
 `uint32_t channel,`
 `array< int32_t >^ amplitude,`
 `array< uint64_t >^ duration,`
 `STG_DestinationEnumNet destType,`
 `bool doClear)`

10.108.2.2 ClearChannelData() `void ClearChannelData (`
 `int channel)`

Delete a Stimulus Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

10.108.2.3 ClearMultiplexedData() `void ClearMultiplexedData ()`

Clears the Stimulation Memory in the STG device.

10.108.2.4 ClearSyncData() `void ClearSyncData (`
 `int channel)`

Delete a Syncout Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

10.108.2.5 CreateSideband() `SidebandData ^ CreateSideband (`
 `array< int32_t >^ StimulusActive,`
 `array< int32_t >^ Syncout,`
 `array< uint64_t >^ Duration,`
 `uint32_t Bit0Time,`
 `uint32_t Bit3Time,`
 `uint32_t Bit4Time)`

Creates the Sideband Channel for the MEA2100 device.

Each datapoint is represented by an signed 32bit integer value. A value 0 means that the stimulation is active during that time. A value 1 means that the stimulation is not active during that time.

The duration is given as a list of 64 bit integers. Durations are given in units of μs . The STG has a resolution of 20 μs .

Parameters

<i>StimulusActive</i>	A list of datapoints which define weather the Stimulus is active or idle at that time as int32.
-----------------------	---

Parameters

<i>Duration</i>	A list of durations as uint64. The time is given in units of μs .
<i>Bit0Time</i>	Time in μs for which Bit 0 (Blanking) is to be extended.

Parameters

<i>Bit3Time</i>	Time in μs for which Bit 3 (Stimulus Enable) is to be extended.
-----------------	--

Parameters

<i>Bit4Time</i>	Time in μs for which Bit 4 (Stimulus Selector) is to be extended.
-----------------	--

Returns

Error Status. 0 on success.

10.108.2.6 ForceStatusEvent() `void ForceStatusEvent ()`

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

10.108.2.7 GetAvailableMemory() `uint32_t GetAvailableMemory ()`

Get the amount of memory available in the currently selected segment of the STG.

Returns

The total memory available on the STG in bytes.

10.108.2.8 GetCurrentRangeInNanoAmp() `int GetCurrentRangeInNanoAmp (
uint32_t channel)`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

10.108.2.9 GetCurrentResolutionInNanoAmp() `int GetCurrentResolutionInNanoAmp (
uint32_t channel)`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

10.108.2.10 GetDACResolution() `int GetDACResolution ()`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

10.108.2.11 GetMultiplexedDataChannelsInBlock() `int GetMultiplexedDataChannelsInBlock ()`

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in `SendMultiplexedData`

10.108.2.12 GetNumberOfAnalogChannels() `int GetNumberOfAnalogChannels ()`

Get the number of STG channels.

Returns

The number of STG channels.

10.108.2.13 GetTotalMemory() `uint32_t GetTotalMemory ()`

Get the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

10.108.2.14 GetVoltageRangeInMicroVolt() `int GetVoltageRangeInMicroVolt (
uint32_t channel)`

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

10.108.2.15 GetVoltageResolutionInMicroVolt() `int GetVoltageResolutionInMicroVolt (
uint32_t channel)`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

10.108.2.16 PrepareAndAppendData() `void PrepareAndAppendData (`
`uint32_t channel,`
`array< int32_t >^ amplitude,`
`array< uint64_t >^ duration,`
`STG_DestinationEnumNet destType)`

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

10.108.2.17 PrepareAndSendData() `void PrepareAndSendData (`
`uint32_t channel,`
`array< int32_t >^ amplitude,`
`array< uint64_t >^ duration,`
`STG_DestinationEnumNet destType)`

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

10.108.2.18 PrepareData() `StimulusDeviceDataAndUnrolledData` ^ PrepareData (
 int *channel*,
 array< int32_t >^ *amplitude*,
 array< uint64_t >^ *duration*,
 STG_DestinationEnumNet *destType*)

10.108.2.19 SendMultiplexedData() void SendMultiplexedData (
 array< uint16_t >^ *data*)

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

Parameters

<i>data</i>	Array of data to be sent.
-------------	---------------------------

10.108.2.20 SendPreparedData() void SendPreparedData (
 int *channel*,
 `StimulusDeviceDataAndUnrolledData`^ *device_data_and_unrolled*,
 STG_DestinationEnumNet *destType*)

10.108.2.21 SendStart() void SendStart (
 uint32_t *triggermap*)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

10.108.2.22 SendStop() [1/2] void SendStop (
 uint32_t *triggermap*)

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

10.108.2.23 SendStop() [2/2] `void SendStop (`
`uint32_t triggermap,`
`int options)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
<i>options</i>	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout associated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active..

10.108.2.24 SetupTrigger() `void SetupTrigger (`
`uint32_t first_trigger,`
`array< uint32_t >^ channelmap,`
`array< uint32_t >^ syncoutmap,`
`array< uint32_t >^ repeat)`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>first_trigger</i>	The number of the first trigger to change.
----------------------	--

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

10.108.2.25 SetupTriggerSingle() `void SetupTriggerSingle (`
`uint32_t trigger,`
`uint32_t channelmap,`
`uint32_t syncoutmap,`
`uint32_t repeat)`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>trigger</i>	The trigger to change.
----------------	------------------------

Parameters

<i>channelmap</i>	A bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	A bitmap of syncouts that belong to this trigger.
<i>repeat</i>	The number of times this trigger should be repeated.

10.108.2.26 StartPoll() `void StartPoll ()`

Starts the interrupt fetching thread and delivers events


10.108.2.27 StopPoll() `void StopPoll ()`

Stops the interrupt fetching thread and delivers events

10.108.3 Event Documentation**10.108.3.1 PollStatusEvent** `OnStgPollStatus^ PollStatusEvent`**10.109 CSw2to64DeviceNet Class Reference**

The class to control the MCS-USB-Sw2to64 device.

Inheritance diagram for CSw2to64DeviceNet:

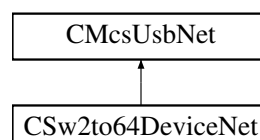

Public Member Functions

- [CSw2to64DeviceNet](#) ()
- [~CSw2to64DeviceNet](#) ()
- unsigned short [GetNumber](#) ()
Gets the number of channels that can be switched in this box.
- array< unsigned char > ^ [GetChannels](#) ()
Gets the current switch positions as char array.
- void [SetChannels](#) (array< unsigned char > ^ pattern)
Sets the switch positions from a char array.
- unsigned char [GetChannel](#) (unsigned short index)
Gets one current switch position.
- void [SetChannel](#) (unsigned short index, unsigned char pattern)
Sets one switch position.

Additional Inherited Members**10.109.1 Detailed Description**

The class to control the MCS-USB-Sw2to64 device.

This class controls the settings of the MCS-USB-Sw2to64. The box has two inputs for signals. Each of the 64 outputs can be connected to one of the input signals, could be held open or connected ground. Valid switch states are 0, 1, 2 or 3 for each of the settings.

10.109.2 Constructor & Destructor Documentation

10.109.2.1 [CSw2to64DeviceNet\(\)](#) [CSw2to64DeviceNet](#) ()

10.109.2.2 [~CSw2to64DeviceNet\(\)](#) [~CSw2to64DeviceNet](#) ()

10.109.3 Member Function Documentation

10.109.3.1 [GetChannel\(\)](#) unsigned char [GetChannel](#) (
 unsigned short *index*)

Gets one current switch position.

Parameters

in	<i>index</i>	number of channel to read the switch position from
----	--------------	--

Returns

switch position of desired channel

10.109.3.2 GetChannels() `array<unsigned char> ^ GetChannels ()`

Gets the current switch positions as char array.

Returns

array of char with the size of the number of channels, each char has the setting of a channel

10.109.3.3 GetNumber() `unsigned short GetNumber ()`

Gets the number of channels that can be switched in this box.

The box can have a different number of channels it can switch. Up to now usually 64 channels are returned

10.109.3.4 SetChannel() `void SetChannel (
 unsigned short index,
 unsigned char pattern)`

Sets one switch position.

Parameters

in	<i>index</i>	number of channel to write the switch position to
in	<i>pattern</i>	switch position of the channel

10.109.3.5 SetChannels() `void SetChannels (
 array< unsigned char >^ pattern)`

Sets the switch positions from a char array.


Parameters

in	<i>pattern</i>	array of char with the size of the number of channels, each char has the setting of a channel
----	----------------	---

10.110 CTcxDeviceNet Class Reference

Class to control a Temperature Controller (TCX)

Inheritance diagram for CTcxDeviceNet:

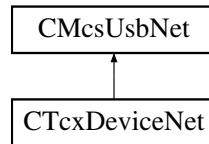

Public Member Functions

- [CTcxDeviceNet](#) ()
Initializes a new instance of [CTcxDeviceNet](#) class.
- [~CTcxDeviceNet](#) ()
- unsigned int [GetNumControlChannels](#) ()
Gets the number of channels the device can control/regulate.
- unsigned int [GetNumMeasureChannels](#) ()
Gets the number of channels the device can measure.
- int [GetValue](#) (unsigned int channel)
Gets the temperate of the specified channel in units of 0.1 °C.
- int [GetValueHires](#) (unsigned int channel)
Gets the temperate of the specified channel in units of 0.01 °C.
- int [GetHeaterTemp](#) (unsigned int channel)
Gets the temperate of the specified heater in units of 0.1 °C.
- int [GetHeaterLimit](#) (unsigned int device)
Gets the temperate limit of the specified heater in units of 0.1 °C.
- double [GetMaxHeaterPowerMultiwell](#) ()
queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W
- void [SetMaxHeaterPowerMultiwell](#) (double MaxPowerWatt)
sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W
- bool [GetHasThermocouple](#) ()
Gets weather the device supports a thermocouple.
- bool [GetEnableHeaterLimit](#) (unsigned int device)
- bool [GetEnableThermocouple](#) (unsigned int device)
- TcxSensorTypeEnumNet [GetSensorType](#) (unsigned int device)
- String ^ [GetUnit](#) (unsigned int channel)
- unsigned int [GetBoardTemp](#) ()
Gets the temperate of the mainboard in units of 0.1 °C.
- unsigned int [GetVoliti](#) (unsigned int channel)
- unsigned int [GetNumDevices](#) ()
- void [SetSetpoint](#) (unsigned int channel, int sp)
Sets the target temperate of specified channel in units of 0.1 °C.
- void [SetDevice](#) (unsigned int channel, int device)
- void [SetOnOff](#) (unsigned int channel, bool on)
Switches the specified channel on or off.
- void [SetCalibration](#) (unsigned int channel, int calib)
- void [SetP](#) (unsigned int device, int p_coeff)
Sets the P-coefficient of the specified device.

- void [SetI](#) (unsigned int device, int i_coeff)
Sets the I-coefficient of the specified device.
- void [SetD](#) (unsigned int device, int d_coeff)
Sets the D-coefficient of the specified device.
- void [SetMaxP](#) (unsigned int device, int maxp)
Sets the maximum heater power of the specified device.
- void [SetHeaterLimit](#) (unsigned int device, int heater_limit)
- void [SetEnableHeaterLimit](#) (unsigned int device, bool enable)
- void [SetEnableThermocouple](#) (unsigned int device, bool enable)
- void [SetSensorType](#) (unsigned int device, TcxSensorTypeEnumNet type)
- void [SetDevname](#) (unsigned int device, String^ Devicename)
- int [GetSetpoint](#) (unsigned int channel)
Gets the target temperate of specified channel in units of 0.1 °C.
- int [GetDevice](#) (unsigned int channel)
- int [GetOnOff](#) (unsigned int channel)
Gets if the specified channel is on or off.
- int [GetCalibration](#) (unsigned int channel)
- int [GetP](#) (unsigned int device)
Gets the P-coefficient of the specified device.
- int [GetI](#) (unsigned int device)
Gets the I-coefficient of the specified device.
- int [GetD](#) (unsigned int device)
Gets the D-coefficient of the specified device.
- int [GetMaxP](#) (unsigned int device)
Gets the maximum heater power of the specified device.
- String ^ [GetDevname](#) (unsigned int device)
- TcxDeviceTypeEnumNet [GetDeviceType](#) ()
- int [GetSetpointMin](#) (unsigned int channel)
- int [GetCalibrationMin](#) (unsigned int channel)
- int [GetPMin](#) (unsigned int device)
- int [GetIMin](#) (unsigned int device)
- int [GetDMin](#) (unsigned int device)
- int [GetMaxpMin](#) (unsigned int device)
- int [GetSetpointMax](#) (unsigned int channel)
- int [GetCalibrationMax](#) (unsigned int channel)
- int [GetPMax](#) (unsigned int device)
- int [GetIMax](#) (unsigned int device)
- int [GetDMax](#) (unsigned int device)
- int [GetMaxpMax](#) (unsigned int device)
- int [GetSetpointDecp](#) (unsigned int channel)
- int [GetCalibrationDecp](#) (unsigned int channel)
- int [GetPDecp](#) (unsigned int device)
- int [GetIDecp](#) (unsigned int device)
- int [GetDDecp](#) (unsigned int device)
- int [GetMaxpDecp](#) (unsigned int device)
- int [GetResX](#) (unsigned int channel)
- int [GetResS](#) (unsigned int channel)
- int [GetRes1](#) (unsigned int channel)
- int [GetRes2](#) (unsigned int channel)
- int [GetPwrSet](#) (unsigned int channel)
- int [GetPwrOut](#) (unsigned int channel)
- int [GetDuty](#) (unsigned int channel)
Gets the duty cycle of the heating element.

- int [GetUOut](#) (unsigned int channel)
Gets the voltage on the heating element.
- int [GetIOut](#) (unsigned int channel)
Gets the current through the heating element.
- int [GetROut](#) (unsigned int channel)
Gets the resistance of the heating element.
- int [GetPOut](#) (unsigned int channel)
Gets the output power of the heating element.
- int [GetCurrent](#) (unsigned int channel)
- int [GetThermocoupleTemp](#) (unsigned int channel)
- int [GetThermocoupleTempAbs](#) (unsigned int channel)
- int [GetThermocoupleReferenceTemp](#) (unsigned int channel)
- unsigned int [GetThermocoupleNanovoltPerKelvin](#) (unsigned int channel)
Gets the proportional constant for the thermocouple.
- void [SetThermocoupleNanovoltPerKelvin](#) (unsigned int channel, unsigned int value)
Sets the proportional constant for the thermocouple.
- int [GetThermocoupleCalibration](#) (unsigned int channel)
- void [CalibrateThermocouple](#) (unsigned int channel)
- void [SetDeviceType](#) (TcxDeviceTypeEnumNet devicetype)
- void [FactoryReset](#) ()

Additional Inherited Members

10.110.1 Detailed Description

Class to control a Temperature Controller (TCX)

10.110.2 Constructor & Destructor Documentation

10.110.2.1 CTcxDeviceNet() [CTcxDeviceNet](#) ()

Initializes a new instance of [CTcxDeviceNet](#) class.

10.110.2.2 ~CTcxDeviceNet() [~CTcxDeviceNet](#) ()

10.110.3 Member Function Documentation

10.110.3.1 CalibrateThermocouple() [void CalibrateThermocouple](#) (unsigned int *channel*)

10.110.3.2 FactoryReset() `void FactoryReset ()`

10.110.3.3 GetBoardTemp() `unsigned int GetBoardTemp ()`

Gets the temperate of the mainboard in units of 0.1 °C.

10.110.3.4 GetCalibration() `int GetCalibration (`
`unsigned int channel)`

10.110.3.5 GetCalibrationDecp() `int GetCalibrationDecp (`
`unsigned int channel)`

10.110.3.6 GetCalibrationMax() `int GetCalibrationMax (`
`unsigned int channel)`

10.110.3.7 GetCalibrationMin() `int GetCalibrationMin (`
`unsigned int channel)`

10.110.3.8 GetCurrent() `int GetCurrent (`
`unsigned int channel)`

10.110.3.9 GetD() `int GetD (`
`unsigned int device)`

Gets the D-coefficient of the specified device.

10.110.3.10 GetDDecp() `int GetDDecp (`
`unsigned int device)`

10.110.3.11 GetDevice() `int GetDevice (`
`unsigned int channel)`

10.110.3.12 GetDeviceType() `TcxDeviceTypeEnumNet GetDeviceType ()`

10.110.3.13 GetDevname() `String ^ GetDevname (`
`unsigned int device)`

10.110.3.14 GetDMax() `int GetDMax (`
`unsigned int device)`

10.110.3.15 GetDMin() `int GetDMin (`
`unsigned int device)`

10.110.3.16 GetDuty() `int GetDuty (`
`unsigned int channel)`

Gets the duty cycle of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The duty cycle in percent, the value of $320 * 64$ corresponds to 100 %.

10.110.3.17 GetEnableHeaterLimit() `bool GetEnableHeaterLimit (`
`unsigned int device)`

10.110.3.18 GetEnableThermocouple() `bool GetEnableThermocouple (`
`unsigned int device)`

10.110.3.19 GetHasThermocouple() `bool GetHasThermocouple ()`

Gets whether the device supports a thermocouple.

10.110.3.20 GetHeaterLimit() `int GetHeaterLimit (
 unsigned int device)`

Gets the temperature limit of the specified heater in units of 0.1 °C.

10.110.3.21 GetHeaterTemp() `int GetHeaterTemp (
 unsigned int channel)`

Gets the temperature of the specified heater in units of 0.1 °C.

10.110.3.22 GetI() `int GetI (
 unsigned int device)`

Gets the I-coefficient of the specified device.

10.110.3.23 GetIDecp() `int GetIDecp (
 unsigned int device)`**10.110.3.24 GetIMax()** `int GetIMax (
 unsigned int device)`**10.110.3.25 GetIMin()** `int GetIMin (
 unsigned int device)`**10.110.3.26 GetIOut()** `int GetIOut (
 unsigned int channel)`

Gets the current through the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The current in units of mA.

10.110.3.27 GetMaxHeaterPowerMultiwell() `double GetMaxHeaterPowerMultiwell ()`

queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

10.110.3.28 GetMaxP() `int GetMaxP (
 unsigned int device)`

Gets the maximum heater power of the specified device.

10.110.3.29 GetMaxpDecp() `int GetMaxpDecp (
 unsigned int device)`

10.110.3.30 GetMaxpMax() `int GetMaxpMax (
 unsigned int device)`

10.110.3.31 GetMaxpMin() `int GetMaxpMin (
 unsigned int device)`

10.110.3.32 GetNumControlChannels() `unsigned int GetNumControlChannels ()`

Gets the number of channels the device can control/regulate.

10.110.3.33 GetNumDevices() `unsigned int GetNumDevices ()`

10.110.3.34 GetNumMeasureChannels() `unsigned int GetNumMeasureChannels ()`

Gets the number of channels the device can measure.

10.110.3.35 GetOnOff() `int GetOnOff (
 unsigned int channel)`

Gets if the specified channel is on or off.

10.110.3.36 GetP() `int GetP (
 unsigned int device)`

Gets the P-coefficient of the specified device.

10.110.3.37 GetPDecp() `int GetPDecp (
 unsigned int device)`**10.110.3.38 GetPMax()** `int GetPMax (
 unsigned int device)`**10.110.3.39 GetPMin()** `int GetPMin (
 unsigned int device)`**10.110.3.40 GetPOut()** `int GetPOut (
 unsigned int channel)`

Gets the output power of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The resistance in units of mW.

10.110.3.41 GetPwrOut() `int GetPwrOut (unsigned int channel)`

10.110.3.42 GetPwrSet() `int GetPwrSet (unsigned int channel)`

10.110.3.43 GetRes1() `int GetRes1 (unsigned int channel)`

10.110.3.44 GetRes2() `int GetRes2 (unsigned int channel)`

10.110.3.45 GetResS() `int GetResS (unsigned int channel)`

10.110.3.46 GetResX() `int GetResX (unsigned int channel)`

10.110.3.47 GetROut() `int GetROut (unsigned int channel)`

Gets the resistance of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The resistance in units of 0.1 Ohm.

10.110.3.48 GetSensorType() `TcxSensorTypeEnumNet GetSensorType (unsigned int device)`

10.110.3.49 GetSetpoint() `int GetSetpoint (unsigned int channel)`

Gets the target temperate of specified channel in units of 0.1 °C.

10.110.3.50 GetSetpointDecp() `int GetSetpointDecp (unsigned int channel)`

10.110.3.51 GetSetpointMax() `int GetSetpointMax (unsigned int channel)`

10.110.3.52 GetSetpointMin() `int GetSetpointMin (unsigned int channel)`

10.110.3.53 GetThermocoupleCalibration() `int GetThermocoupleCalibration (unsigned int channel)`

10.110.3.54 GetThermocoupleNanovoltPerKelvin() `unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)`

Gets the proportional constant for the thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

10.110.3.55 GetThermocoupleReferenceTemp() `int GetThermocoupleReferenceTemp (unsigned int channel)`

10.110.3.56 GetThermocoupleTemp() `int GetThermocoupleTemp (unsigned int channel)`

10.110.3.57 GetThermocoupleTempAbs() `int GetThermocoupleTempAbs (unsigned int channel)`

10.110.3.58 GetUnit() `String ^ GetUnit (unsigned int channel)`

10.110.3.59 GetUOut() `int GetUOut (unsigned int channel)`

Gets the voltage on the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The voltage in units of mV.

10.110.3.60 GetValue() `int GetValue (unsigned int channel)`

Gets the temperate of the specified channel in units of 0.1 °C.

10.110.3.61 GetValueHires() `int GetValueHires (unsigned int channel)`

Gets the temperate of the specified channel in units of 0.01 °C.

10.110.3.62 GetVolti() `unsigned int GetVolti (unsigned int channel)`

10.110.3.63 SetCalibration() void SetCalibration (
 unsigned int *channel*,
 int *calib*)

10.110.3.64 SetD() void SetD (
 unsigned int *device*,
 int *d_coeff*)

Sets the D-coefficient of the specified device.

10.110.3.65 SetDevice() void SetDevice (
 unsigned int *channel*,
 int *device*)

10.110.3.66 SetDeviceType() void SetDeviceType (
 TcxDeviceTypeEnumNet *devicetype*)

10.110.3.67 SetDevname() void SetDevname (
 unsigned int *device*,
 String^ *Devicename*)

10.110.3.68 SetEnableHeaterLimit() void SetEnableHeaterLimit (
 unsigned int *device*,
 bool *enable*)

10.110.3.69 SetEnableThermocouple() void SetEnableThermocouple (
 unsigned int *device*,
 bool *enable*)

10.110.3.70 SetHeaterLimit() void SetHeaterLimit (
 unsigned int *device*,
 int *heater_limit*)

10.110.3.71 SetI() void SetI (
 unsigned int *device*,
 int *i_coeff*)

Sets the I-coefficient of the specified device.

10.110.3.72 SetMaxHeaterPowerMultiwell() void SetMaxHeaterPowerMultiwell (
 double *MaxPowerWatt*)

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

10.110.3.73 SetMaxP() void SetMaxP (
 unsigned int *device*,
 int *maxp*)

Sets the maximum heater power of the specified device.

10.110.3.74 SetOnOff() void SetOnOff (
 unsigned int *channel*,
 bool *on*)

Switches the specified channel on or off.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

10.110.3.75 SetP() void SetP (
 unsigned int *device*,
 int *p_coeff*)

Sets the P-coefficient of the specified device.

10.110.3.76 SetSensorType() void SetSensorType (
 unsigned int *device*,
 TcxSensorTypeEnumNet *type*)

10.110.3.77 SetSetpoint() `void SetSetpoint (`
 `unsigned int channel,`
 `int sp)`

Sets the target temperate of specified channel in units of 0.1 °C.

10.110.3.78 SetThermocoupleNanovoltPerKelvin() `void SetThermocoupleNanovoltPerKelvin (`
 `unsigned int channel,`
 `unsigned int value)`

Sets the proportional constant for the thermocouple.


Parameters

<i>channel</i>	Thermocouple channel number.
<i>value</i>	Proportional constant in Nanovolt per Kelvin.

10.111 CTEERFunctionNet Class Reference

[CTEERFunctionNet](#) is the class to control the TEER device

Inheritance diagram for CTEERFunctionNet:

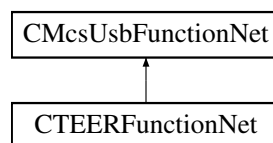

Public Member Functions

- [CTEERFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pTEERFunctionPointer↔ Container)
Initializes a new instance of the [CTEERFunctionNet](#) class.
- [CTEERFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CTEERFunctionNet](#) ()
- [!CTEERFunctionNet](#) ()
- `uint32_t` [GetPeriod_us](#) ()
gets the period of TEER stimulation in us
- void [SetPeriod_us](#) (`uint32_t` period_us)
sets the period of TEER stimulation in us
- `uint32_t` [GetAmplitude_nA](#) ()
gets TEER stimulation amplitude in nA
- void [SetAmplitude_nA](#) (`uint32_t` Amplitude_nA)
sets TEER stimulation amplitude in nA
- TeerWaveformEnumNet [GetWaveform](#) ()
gets TEER stimulation waveform (sine/rect)

- void [SetWaveform](#) (TeerWaveformEnumNet Waveform)
sets TEER stimulation waveform (sine/rect)
- TeerClampModeEnumNet [GetClampMode](#) ()
gets TEER clamp mode (voltage/current)
- void [SetClampMode](#) (TeerClampModeEnumNet ClampMode)
sets TEER clamp mode (voltage/current)
- void [StartSampling](#) (uint32_t NumberOfCycles)
starts TEER stimulation (duration: n cycles) and samples during last cycle
- void [StopSampling](#) ()
stops TEER stimulation and sampling
- uint32_t [IsSamplingFinished](#) ()
returns false iff stimulation/sampling is going on, otherwise true
- void [SetControllerParams](#) (uint32_t P, uint32_t I, uint32_t D)
sets PID controller parameters for voltage clamp mode
- void [GetControllerParams](#) ([System::Runtime::InteropServices::Out]uint32_t% P, [System::Runtime::InteropServices::Out]uint32_t% I, [System::Runtime::InteropServices::Out]uint32_t% D)
gets PID controller parameters for voltage clamp mode
- array< int32_t > ^ [GetSampleBufferChunk](#) (int Buffer_Length)
private function to query max. 100 bytes of sample buffer; called internally
- array< int32_t > ^ [GetSampleVoltageBuffer_uV](#) (int Buffer_Length)
returns voltage sample buffer (max. 500 values); unit: uV
- uint32_t [GetMaxChunkSize_Byte](#) ()
private function to be called internally only
- uint32_t [GetBytesPerSample](#) ()
private function to be called internally only
- uint32_t [GetNumberOfAvailableSamples](#) ()
private function to be called internally only
- void [SetBufferIndex](#) (uint32_t NewBufferIndex)
pre-selects sample buffer to be tranferred by [GetSampleVoltageBuffer_uV\(\)](#)
- uint32_t [GetAdapterCode](#) ()
- uint32_t [GetRotaryPositionCode](#) ()
- void [SetExternalLED](#) (uint32_t NewState)
- void [SetCurrentEnable](#) (bool NewCurrentEnable)
when disabled, no current will flow through chamber
- bool [GetCurrentEnable](#) ()
when disabled, no current will flow through chamber
- int32_t [GetUptimeSeconds](#) ()
returns time in seconds since device was powered up
- void [StartInternalCalibration](#) ()
starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call
- bool [IsInternalCalibrationFinished](#) ()
queries whether internal calibration has finished
- int [GetDacZero](#) ()
returns DAC-offset (result of internal calibration); use to check for plausibility only
- void [CancelInternalCalibration](#) ()
in case the internal calibration "hangs", this will cancel it
- void [SetLiquidResistance](#) (int32_t NewLiquidResistance_Ohm)
- int32_t [GetLiquidResistance](#) ()
- int [GetScaleFactorU1](#) ()
returns U1 scale factor times 10⁶ (result of internal calibration)
- int [GetScaleFactorU2](#) ()

- returns U2 scale factor times 10^6 (result of internal calibration)*
 • int [GetAdcOffsetU1](#) ()
returns ADC offset of U1 channel (result of internal calibration)
- int [GetAdcOffsetU2](#) ()
returns ADC offset of U2 channel (result of internal calibration)

Additional Inherited Members

10.111.1 Detailed Description

[CTEERFunctionNet](#) is the class to control the TEER device

10.111.2 Constructor & Destructor Documentation

10.111.2.1 CTEERFunctionNet() [1/2] [CTEERFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pTEERFunctionPointerContainer*)

Initializes a new instance of the [CTEERFunctionNet](#) class.

10.111.2.2 CTEERFunctionNet() [2/2] [CTEERFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

10.111.2.3 ~CTEERFunctionNet() [virtual ~CTEERFunctionNet](#) () [virtual]

10.111.2.4 "!CTEERFunctionNet() [!CTEERFunctionNet](#) ()

10.111.3 Member Function Documentation

10.111.3.1 CancelInternalCalibration() [void CancelInternalCalibration](#) ()

in case the internal calibration "hangs", this will cancel it

10.111.3.2 GetAdapterCode() `uint32_t GetAdapterCode ()`

Returns

10.111.3.3 GetAdcOffsetU1() `int GetAdcOffsetU1 ()`

returns ADC offset of U1 channel (result of internal calibration)

Returns

10.111.3.4 GetAdcOffsetU2() `int GetAdcOffsetU2 ()`

returns ADC offset of U2 channel (result of internal calibration)

Returns

10.111.3.5 GetAmplitude_nA() `uint32_t GetAmplitude_nA ()`

gets TEER stimulation amplitude in nA

Returns

current stimulation amplitude in nA

10.111.3.6 GetBytesPerSample() `uint32_t GetBytesPerSample ()`

private function to be called internally only

Returns

10.111.3.7 GetClampMode() `TeerClampModeEnumNet GetClampMode ()`

gets TEER clamp mode (voltage/current)

Returns

current TEER clamp mode

10.111.3.8 GetControllerParams() `void GetControllerParams ([System::Runtime::InteropServices::Out] uint32_t% P, [System::Runtime::InteropServices::Out] uint32_t% I, [System::Runtime::InteropServices::Out] uint32_t% D)`

gets PID controller parameters for voltage clamp mode

Parameters

<i>P</i>	
<i>I</i>	
<i>D</i>	

10.111.3.9 GetCurrentEnable() `bool GetCurrentEnable ()`

when disabled, no current will flow through chamber

Returns

10.111.3.10 GetDacZero() `int GetDacZero ()`

returns DAC-offset (result of internal calibration); use to check for plausibility only

Returns

10.111.3.11 GetLiquidResistance() `int32_t GetLiquidResistance ()`

Returns

10.111.3.12 GetMaxChunkSize_Byte() `uint32_t GetMaxChunkSize_Byte ()`

private function to be called internally only

Returns

10.111.3.13 GetNumberOfAvailableSamples() `uint32_t GetNumberOfAvailableSamples ()`

private function to be called internally only

Returns

10.111.3.14 GetPeriod_us() `uint32_t GetPeriod_us ()`

gets the period of TEER stimulation in us

Returns

10.111.3.15 GetRotaryPositionCode() `uint32_t GetRotaryPositionCode ()`

Returns

10.111.3.16 GetSampleBufferChunk() `array<int32_t> ^ GetSampleBufferChunk (
int Buffer_Length)`

private function to query max. 100 bytes of sample buffer; called internally

Parameters

<i>Buffer_Length</i>	The maximal length of Buffer.
----------------------	-------------------------------

Returns

10.111.3.17 GetSampleVoltageBuffer_uV() `array<int32_t> ^ GetSampleVoltageBuffer_uV (
int Buffer_Length)`

returns voltage sample buffer (max. 500 values); unit: uV

Parameters

<i>Buffer_Length</i>	The maximal length of Buffer.
----------------------	-------------------------------

Returns

10.111.3.18 GetScaleFactorU1() `int GetScaleFactorU1 ()`

returns U1 scale factor times 10^6 (result of internal calibration)

Returns

10.111.3.19 GetScaleFactorU2() `int GetScaleFactorU2 ()`

returns U2 scale factor times 10^6 (result of internal calibration)

Returns

10.111.3.20 GetUptimeSeconds() `int32_t GetUptimeSeconds ()`

returns time in seconds since device was powered up

Returns

10.111.3.21 GetWaveform() `TeerWaveformEnumNet GetWaveform ()`

gets TEER stimulation waveform (sine/rect)

Returns

waveform enum

10.111.3.22 IsInternalCalibrationFinished() `bool IsInternalCalibrationFinished ()`

queries whether internal calibration has finished

Returns

10.111.3.23 IsSamplingFinished() `uint32_t IsSamplingFinished ()`

returns false iff stimulation/sampling is going on, otherwise true

Returns

10.111.3.24 SetAmplitude_nA() `void SetAmplitude_nA (
uint32_t Amplitude_nA)`

sets TEER stimulation amplitude in nA

Parameters

<i>Amplitude_nA</i>	new stimulation amplitude in nA
---------------------	---------------------------------

10.111.3.25 SetBufferIndex() `void SetBufferIndex (
uint32_t NewBufferIndex)`

pre-selects sample buffer to be tranferred by [GetSampleVoltageBuffer_uV\(\)](#)

Parameters

<i>NewBufferIndex</i>	0 - chamber voltage; 1 - compliance voltage
-----------------------	---

10.111.3.26 SetClampMode() `void SetClampMode (
TeerClampModeEnumNet ClampMode)`

sets TEER clamp mode (voltage/current)

Parameters

<i>ClampMode</i>	new TEER clamp mode
------------------	---------------------

10.111.3.27 SetControllerParams() void SetControllerParams (
 uint32_t *P*,
 uint32_t *I*,
 uint32_t *D*)

sets PID controller parameters for voltage clamp mode

Parameters

<i>P</i>	
<i>I</i>	
<i>D</i>	

10.111.3.28 SetCurrentEnable() void SetCurrentEnable (
 bool *NewCurrentEnable*)

when disabled, no current will flow through chamber

Parameters

<i>NewCurrentEnable</i>	
-------------------------	--

10.111.3.29 SetExternalLED() void SetExternalLED (
 uint32_t *NewState*)

Parameters

<i>NewState</i>	
-----------------	--

10.111.3.30 SetLiquidResistance() void SetLiquidResistance (
 int32_t *NewLiquidResistance_Ohm*)

Parameters

<i>NewLiquidResistance_Ohm</i>	
--------------------------------	--

10.111.3.31 SetPeriod_us() `void SetPeriod_us (`
`uint32_t period_us)`

sets the period of TEER stimulation in us

Parameters

<i>period_us</i>	
------------------	--

10.111.3.32 SetWaveform() `void SetWaveform (`
`TeerWaveformEnumNet Waveform)`

sets TEER stimulation waveform (sine/rect)

Parameters

<i>Waveform</i>	waveform enum
-----------------	---------------

10.111.3.33 StartInternalCalibration() `void StartInternalCalibration ()`

starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call

10.111.3.34 StartSampling() `void StartSampling (`
`uint32_t NumberOfCycles)`

starts TEER stimulation (duration: n cycles) and samples during last cycle

Parameters


<i>NumberOfCycles</i>	number of cycles (sine or rect) to output (0 - loop forever)
-----------------------	--

10.111.3.35 StopSampling() `void StopSampling ()`

stops TEER stimulation and sampling

10.112 CTEERMachineDeviceNet Class Reference

Inheritance diagram for CTEERMachineDeviceNet:

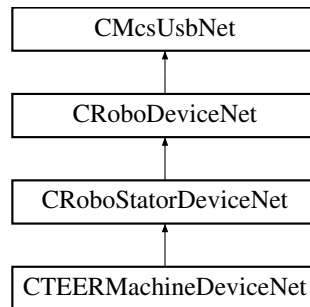

Public Member Functions

- [CTEERMachineDeviceNet](#) ()
- [~CTEERMachineDeviceNet](#) ()

Properties

- [CTEERFunctionNet](#)[^] [TEERFunctionNet](#) [get]

Additional Inherited Members

10.112.1 Constructor & Destructor Documentation

10.112.1.1 [CTEERMachineDeviceNet](#)() [CTEERMachineDeviceNet](#) ()

10.112.1.2 [~CTEERMachineDeviceNet](#)() [~CTEERMachineDeviceNet](#) ()


10.112.2 Property Documentation

10.112.2.1 [TEERFunctionNet](#) [CTEERFunctionNet](#)[^] [TEERFunctionNet](#) [get]

10.113 CUsbExceptionNet Class Reference

Exception class that is thrown in case of an USB error.

Inheritance diagram for CUsbExceptionNet:

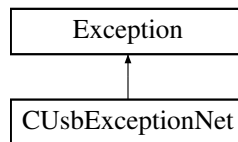

Public Member Functions

- [CUsbExceptionNet](#) (uint32_t status)
Constructor of a CUsbException.
- [CUsbExceptionNet](#) (uint32_t status, String^ message)

Properties

- uint32_t [Status](#) [get]

10.113.1 Detailed Description

Exception class that is thrown in case of an USB error.

10.113.2 Constructor & Destructor Documentation

10.113.2.1 CUsbExceptionNet() [1/2] [CUsbExceptionNet](#) (uint32_t status)

Constructor of a CUsbException.

Parameters

<i>status</i>	the status number
---------------	-------------------

10.113.2.2 CUsbExceptionNet() [2/2] [CUsbExceptionNet](#) (uint32_t status, String^ message)

10.113.3 Property Documentation

10.113.3.1 Status `uint32_t Status [get]`

10.114 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference

Public Member Functions

- [CVoltageRangeInfoNet](#) (int vr, String^ vrString)

Public Attributes

- int [VoltageRangeInMicroVolt](#)
- String ^ [VoltageRangeDisplayStringMilliVolt](#)

10.114.1 Constructor & Destructor Documentation

10.114.1.1 CVoltageRangeInfoNet() `CVoltageRangeInfoNet (`
 int vr,
 String^ vrString)


10.114.2 Member Data Documentation

10.114.2.1 VoltageRangeDisplayStringMilliVolt `String ^ VoltageRangeDisplayStringMilliVolt`

10.114.2.2 VoltageRangeInMicroVolt `int VoltageRangeInMicroVolt`

10.115 CW2100_FunctionNet Class Reference

Inheritance diagram for CW2100_FunctionNet:

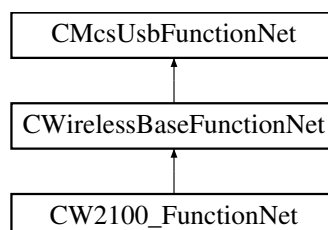

Classes

- struct [AudioChannelsNet](#)

Public Member Functions

- [CW2100_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] w2100_FunctionPointerContainer)
- [CW2100_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- [array< HeadStageIDType[^]> ^ GetAvailableHeadstages](#) (unsigned int max_length)
- [void SelectHeadstage](#) (unsigned int IDorEntry, int TimeSlotNr)
- [void DeselectHeadstage](#) (int TimeSlotNr)
- [void DeselectAllHeadstages](#) ()
- [HeadStageIDTypeState[^] GetSelectedHeadstageState](#) (int TimeSlotNr)
- [BatteryState[^] GetBatteryState](#) (int TimeSlotNr)
- [System::String[^] GetUserDefinedName](#) (unsigned short ID)
- [System::String[^] GetUserDefinedNameFromSelectedHS](#) (int TimeSlotNr)
- [System::String[^] GetUserDefinedNameCache](#) (unsigned short ID)
- [uint32_t GetUserDefinedNameCache](#) (unsigned short ID, [System::Runtime::InteropServices::Out]System::String[^]% Name)
- [W2100_StimulusParametersNet[^] GetStimulusParameters](#) (unsigned short ID)
- [W2100_StimulusParametersNet[^] GetStimulusParametersFromSelectedHS](#) (int TimeSlotNr)
- [W2100_StimulusParametersNet[^] GetStimulusParametersCache](#) (unsigned int typeValue)
- [uint32_t GetStimulusParametersCache](#) (unsigned int typeValue, [System::Runtime::InteropServices::Out]W2100_StimulusParametersNet[^]% StimulusParameters)
- [void SetSelectedChannels](#) (array< BYTE >[^] channels, int TimeSlotNr)
- [array< BYTE > ^ GetSelectedChannels](#) (int TimeSlotNr)
- [void SetMultiHeadstageMode](#) (bool Mode)
- [bool GetMultiHeadstageMode](#) ()
- [void SetHeadstageSamplingActive](#) (bool Active, int TimeSlotNr)
- [bool GetHeadstageSamplingActive](#) (int TimeSlotNr)
- [void SetHeadstageToSleep](#) (unsigned int Sleep16ms, int TimeSlotNr)
- [void SetHeadstageOnOff](#) (unsigned short On, int TimeSlotNr)
- [unsigned short GetHeadstageOnOff](#) (int TimeSlotNr)
- [unsigned int GetAnalogOutChannel](#) ([System::Runtime::InteropServices::Out]int % automatic, unsigned short index)
- [void SetAnalogOutChannel](#) (int automatic, unsigned short index, unsigned int Channel)
- [array< unsigned int > ^ GetAnalogOutFilter](#) ([System::Runtime::InteropServices::Out]int % automatic)
- [void SetAnalogOutFilter](#) (int automatic, array< unsigned int >[^] Coeffs)
- [AnalogOut_DAC_Range_EnumNet GetDacRange](#) ()
- [void SetDacRange](#) (AnalogOut_DAC_Range_EnumNet range)
- [CFilterPropertyNet[^] GetFilterProperty](#) (W2100DacqGroupChannelEnumNet GroupID, unsigned int index)
- [array< CFilterPropertyNet[^]> ^ GetFilterProperties](#) (W2100DacqGroupChannelEnumNet GroupID)
- [void SetAccelGyroEnabled](#) (W2100_Accel_Gyro_Select_EnumNet enable, int TimeSlotNr)
- [W2100_Accel_Gyro_Select_EnumNet GetAccelGyroEnabled](#) (int TimeSlotNr)
- [void SetAccelGyroDesiredRate](#) (int rate, int TimeSlotNr)
- [int GetAccelGyroDesiredRate](#) (int TimeSlotNr)
- [int GetAccelGyroCurrentRate](#) (int TimeSlotNr)
- [void SetAccelRange](#) (int range, int TimeSlotNr)
- [int GetAccelRange](#) (int TimeSlotNr)
- [void SetGyroRange](#) (int range, int TimeSlotNr)
- [int GetGyroRange](#) (int TimeSlotNr)
- [void SetAudioChannels](#) (array< [AudioChannelsNet](#)[^]>[^] channels)
- [array< AudioChannelsNet[^]> ^ GetAudioChannels](#) ()
- [unsigned int GetPicFirmwareType](#) (int TimeSlotNr)
- [unsigned int GetFPGA FirmwareType](#) (int TimeSlotNr)

Static Public Member Functions

- static void [ClearUserDefinedNameCache](#) ()
- static void [ClearUserDefinedNameCache](#) (unsigned short ID)
- static void [ClearStimulusParametersCache](#) ()
- static void [ClearStimulusParametersCache](#) (unsigned short ID)

Properties

- [CW2100_StimulatorFunctionNet^ Stimulator](#) [get]
- [CPulseGeneratorFunctionNet^ PulseGenerator](#) [get]

Additional Inherited Members

10.115.1 Constructor & Destructor Documentation

10.115.1.1 CW2100_FunctionNet() [1/2] [CW2100_FunctionNet](#) (
[CMcsUsbNet^ mcsusb](#),
[CMcsUsbFunctionPointerContainer^ w2100_FunctionPointerContainer](#))

10.115.1.2 CW2100_FunctionNet() [2/2] [CW2100_FunctionNet](#) (
[CMcsUsbNet^ mcsusb](#))

10.115.2 Member Function Documentation

10.115.2.1 ClearStimulusParametersCache() [1/2] static void [ClearStimulusParametersCache](#) ()
[static]

10.115.2.2 ClearStimulusParametersCache() [2/2] static void [ClearStimulusParametersCache](#) (
unsigned short *ID*) [static]

10.115.2.3 ClearUserDefinedNameCache() [1/2] static void [ClearUserDefinedNameCache](#) () [static]

10.115.2.4 ClearUserDefinedNameCache() [2/2] static void ClearUserDefinedNameCache (unsigned short *ID*) [static]

10.115.2.5 DeselectAllHeadstages() void DeselectAllHeadstages ()

10.115.2.6 DeselectHeadstage() void DeselectHeadstage (int *TimeSlotNr*)

10.115.2.7 GetAccelGyroCurrentRate() int GetAccelGyroCurrentRate (int *TimeSlotNr*)

10.115.2.8 GetAccelGyroDesiredRate() int GetAccelGyroDesiredRate (int *TimeSlotNr*)

10.115.2.9 GetAccelGyroEnabled() W2100_Accel_Gyro_Select_EnumNet GetAccelGyroEnabled (int *TimeSlotNr*)

10.115.2.10 GetAccelRange() int GetAccelRange (int *TimeSlotNr*)

10.115.2.11 GetAnalogOutChannel() unsigned int GetAnalogOutChannel ([System::Runtime::InteropServices::Out] int % *automatic*, unsigned short *index*)

10.115.2.12 GetAnalogOutFilter() array<unsigned int> ^ GetAnalogOutFilter ([System::Runtime::InteropServices::Out] int % *automatic*)

10.115.2.13 GetAudioChannels() array<AudioChannelsNet^> ^ GetAudioChannels ()

10.115.2.14 GetAvailableHeadstages() array<HeadStageIDType^> ^ GetAvailableHeadstages (
 unsigned int *max_length*)

10.115.2.15 GetBatteryState() BatteryState ^ GetBatteryState (
 int *TimeSlotNr*)

10.115.2.16 GetDacRange() AnalogOut_DAC_Range_EnumNet GetDacRange ()

10.115.2.17 GetFilterProperties() array<CFilterPropertyNet^> ^ GetFilterProperties (
 W2100DacqGroupChannelEnumNet *GroupID*)

10.115.2.18 GetFilterProperty() CFilterPropertyNet ^ GetFilterProperty (
 W2100DacqGroupChannelEnumNet *GroupID*,
 unsigned int *index*)

10.115.2.19 GetFPGA FirmwareType() unsigned int GetFPGA FirmwareType (
 int *TimeSlotNr*)

10.115.2.20 GetGyroRange() int GetGyroRange (
 int *TimeSlotNr*)

10.115.2.21 GetHeadstageOnOff() unsigned short GetHeadstageOnOff (
 int *TimeSlotNr*)

10.115.2.22 GetHeadstageSamplingActive() bool GetHeadstageSamplingActive (
 int *TimeSlotNr*)

10.115.2.23 GetMultiHeadstageMode() bool GetMultiHeadstageMode ()

10.115.2.24 GetPicFirmwareType() unsigned int GetPicFirmwareType (
int *TimeSlotNr*)

10.115.2.25 GetSelectedChannels() array<BYTE> ^ GetSelectedChannels (
int *TimeSlotNr*)

10.115.2.26 GetSelectedHeadstageState() HeadStageIDTypeState ^ GetSelectedHeadstageState (
int *TimeSlotNr*)

10.115.2.27 GetStimulusParametersCache() [1/2] W2100_StimulusParametersNet ^ GetStimulus↔
ParametersCache (
unsigned int *typeValue*)

10.115.2.28 GetStimulusParametersCache() [2/2] uint32_t GetStimulusParametersCache (
unsigned int *typeValue*,
[System::Runtime::InteropServices::Out] W2100_StimulusParametersNet ^% *Stimulus↔*
Parameters)

10.115.2.29 GetStimulusParametersFromSelectedHS() W2100_StimulusParametersNet ^ GetStimulus↔
ParametersFromSelectedHS (
int *TimeSlotNr*)

10.115.2.30 GetStiumlusParameters() W2100_StimulusParametersNet ^ GetStiumlusParameters (
unsigned short *ID*)

10.115.2.31 GetUserDefinedName() System::String ^ GetUserDefinedName (
unsigned short *ID*)

10.115.2.32 GetUserDefinedNameCache() [1/2] System::String ^ GetUserDefinedNameCache (
unsigned short *ID*)

10.115.2.33 GetUserDefinedNameCache() [2/2] uint32_t GetUserDefinedNameCache (unsigned short *ID*, [System::Runtime::InteropServices::Out] System::String^% *Name*)

10.115.2.34 GetUserDefinedNameFromSelectedHS() System::String ^ GetUserDefinedNameFromSelectedHS (int *TimeSlotNr*)

10.115.2.35 SelectHeadstage() void SelectHeadstage (unsigned int *IDorEntry*, int *TimeSlotNr*)

10.115.2.36 SetAccelGyroDesiredRate() void SetAccelGyroDesiredRate (int *rate*, int *TimeSlotNr*)

10.115.2.37 SetAccelGyroEnabled() void SetAccelGyroEnabled (W2100_Accel_Gyro_Select_EnumNet *enable*, int *TimeSlotNr*)

10.115.2.38 SetAccelRange() void SetAccelRange (int *range*, int *TimeSlotNr*)

10.115.2.39 SetAnalogOutChannel() void SetAnalogOutChannel (int *automatic*, unsigned short *index*, unsigned int *Channel*)

10.115.2.40 SetAnalogOutFilter() void SetAnalogOutFilter (int *automatic*, array< unsigned int >^ *Coeffs*)

10.115.2.41 SetAudioChannels() void SetAudioChannels (
array< AudioChannelsNet^>^ channels)

10.115.2.42 SetDacRange() void SetDacRange (
AnalogOut_DAC_Range_EnumNet range)

10.115.2.43 SetGyroRange() void SetGyroRange (
int range,
int TimeSlotNr)

10.115.2.44 SetHeadstageOnOff() void SetHeadstageOnOff (
unsigned short On,
int TimeSlotNr)

10.115.2.45 SetHeadstageSamplingActive() void SetHeadstageSamplingActive (
bool Active,
int TimeSlotNr)

10.115.2.46 SetHeadstageToSleep() void SetHeadstageToSleep (
unsigned int Sleep16ms,
int TimeSlotNr)

10.115.2.47 SetMultiHeadstageMode() void SetMultiHeadstageMode (
bool Mode)

10.115.2.48 SetSelectedChannels() void SetSelectedChannels (
array< BYTE >^ channels,
int TimeSlotNr)


10.115.3 Property Documentation

10.115.3.1 PulseGenerator [CPulseGeneratorFunctionNet](#)[^] PulseGenerator [get]

10.115.3.2 Stimulator [CW2100_StimulatorFunctionNet](#)[^] Stimulator [get]

10.116 CW2100_StimulatorFunctionNet Class Reference

Inheritance diagram for CW2100_StimulatorFunctionNet:

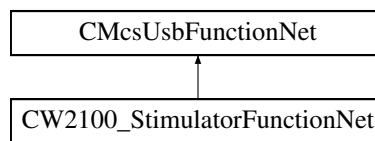

Public Member Functions

- [CW2100_StimulatorFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] [PrepareData](#) (int channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, STG_DestinationEnumNet destType, uint32_t repeat)
- [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] [PrepareDataSync](#) (int channel, array< int32_t >[^] amplitude, array< uint32_t >[^] Sync, array< uint64_t >[^] duration, STG_DestinationEnumNet destType, uint32_t repeat)
- void [SendPreparedData](#) (int channel, [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] device_data_and_unrolled, STG_DestinationEnumNet destType)
- void [ClearChannelData](#) (int channel)
Delete a Stimulus Pattern from STG memory
- int [GetDACResolution](#) ()
Gets number of bits of the DAC resolution.
- int [GetTimeResolutionInNanoSeconds](#) ()
Gets number of bits of the DAC resolution.
- int [GetVoltageRangeInMicroVolt](#) (uint32_t channel)
Gets the Voltage Range of the specified channel in Microvolts.
- int [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)
Gets the Voltage Resolution of the specified channel in Microvolts.
- int [GetCurrentRangeInNanoAmp](#) (uint32_t channel)
Gets the Current Range of the specified channel in Nanoamps.
- int [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)
Gets the Current Resolution of the specified channel in Nanoamps.
- uint32_t [GetNumberOfAnalogChannels](#) ()
- uint32_t [GetNumberOfSyncoutChannels](#) ()
- uint32_t [GetNumberOfTriggerInputs](#) ()
- void [SelectTimeSlot](#) (int TimeSlotNr)
- int [GetTimeSlot](#) ()
- uint32_t [GetStimulationPatternMemory](#) ()

- uint32_t [GetBoostPreTime](#) ()
- void [SetBoostAlwaysOnMode](#) (uint32_t alwayson_mode)
- uint32_t [GetBoostAlwaysOnMode](#) ()
- void [SetDigitalStimulatorTrigger](#) (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, W2100DigitalSourceEnumNet digstream_source, int bitnumber_offset)
- void [GetDigitalStimulatorTrigger](#) (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% digstream_source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
- void [SetDigitalStimulatorTriggerSlope](#) (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, DigitalStimulatorTriggerSlopeEnumNet slope)
- DigitalStimulatorTriggerSlopeEnumNet [GetDigitalStimulatorTriggerSlope](#) (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number)
- void [StartPoll](#) ()
- void [StopPoll](#) ()

Static Public Attributes

- static const uint32_t [BOOST_BIT](#) = (1 << 0)
- static const uint32_t [GND_SWITCH_BIT](#) = (1 << 1)
- static const uint32_t [SYNC_BIT0](#) = (1 << 2)
- static const uint32_t [SYNC_BIT1](#) = (1 << 3)

Events

- [OnStgPollStatus](#)[^] [PollStatusEvent](#)

Additional Inherited Members

10.116.1 Constructor & Destructor Documentation

10.116.1.1 CW2100_StimulatorFunctionNet() [CW2100_StimulatorFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

10.116.2 Member Function Documentation

10.116.2.1 ClearChannelData() void [ClearChannelData](#) (
int *channel*)

Delete a Stimulus Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

10.116.2.2 GetBoostAlwaysOnMode() `uint32_t GetBoostAlwaysOnMode ()`

10.116.2.3 GetBoostPreTime() `uint32_t GetBoostPreTime ()`

10.116.2.4 GetCurrentRangeInNanoAmp() `int GetCurrentRangeInNanoAmp (
uint32_t channel)`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

10.116.2.5 GetCurrentResolutionInNanoAmp() `int GetCurrentResolutionInNanoAmp (
uint32_t channel)`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

10.116.2.6 GetDACResolution() `int GetDACResolution ()`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

10.116.2.7 GetDigitalStimulatorTrigger() void GetDigitalStimulatorTrigger (
 int *TimeSlotNr*,
 DigitalStimulatorTriggerEventEnumNet *trigger_event*,
 int *trigger_number*,
 [System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% *digstream_↵*
source,
 [System::Runtime::InteropServices::Out] int% *bitnumber_offset*)

10.116.2.8 GetDigitalStimulatorTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetDigital↵
 StimulatorTriggerSlope (
 int *TimeSlotNr*,
 DigitalStimulatorTriggerEventEnumNet *trigger_event*,
 int *trigger_number*)

10.116.2.9 GetNumberOfAnalogChannels() uint32_t GetNumberOfAnalogChannels ()

10.116.2.10 GetNumberOfSyncoutChannels() uint32_t GetNumberOfSyncoutChannels ()

10.116.2.11 GetNumberOfTriggerInputs() uint32_t GetNumberOfTriggerInputs ()

10.116.2.12 GetStimulationPatternMemory() uint32_t GetStimulationPatternMemory ()

10.116.2.13 GetTimeResolutionInNanoSeconds() int GetTimeResolutionInNanoSeconds ()

Gets number of bits of the DAC resolution.

Returns

The time resolution in ns.

10.116.2.14 GetTimeSlot() int GetTimeSlot ()

10.116.2.15 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (
 uint32_t *channel*)

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

10.116.2.16 GetVoltageResolutionInMicroVolt() `int GetVoltageResolutionInMicroVolt (uint32_t channel)`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

10.116.2.17 PrepareData() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareData (`

```
int channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType,
uint32_t repeat )
```

10.116.2.18 PrepareDataSync() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareDataSync (`

```
int channel,
array< int32_t >^ amplitude,
array< uint32_t >^ Sync,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType,
uint32_t repeat )
```

10.116.2.19 SelectTimeSlot() `void SelectTimeSlot (int TimeSlotNr)`

10.116.2.20 SendPreparedData() void SendPreparedData (
 int channel,
 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled,
 STG_DestinationEnumNet destType)

10.116.2.21 SendStart() void SendStart (
 uint32_t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

10.116.2.22 SendStop() void SendStop (
 uint32_t triggermap)

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

10.116.2.23 SetBoostAlwaysOnMode() void SetBoostAlwaysOnMode (
 uint32_t alwayson_mode)

10.116.2.24 SetDigitalStimulatorTrigger() void SetDigitalStimulatorTrigger (
 int TimeSlotNr,
 DigitalStimulatorTriggerEventEnumNet trigger_event,
 int trigger_number,
 W2100DigitalSourceEnumNet digstream_source,
 int bitnumber_offset)

10.116.2.25 SetDigitalStimulatorTriggerSlope() void SetDigitalStimulatorTriggerSlope (
 int TimeSlotNr,
 DigitalStimulatorTriggerEventEnumNet trigger_event,
 int trigger_number,
 DigitalStimulatorTriggerSlopeEnumNet slope)

10.116.2.26 StartPoll() `void StartPoll ()`

10.116.2.27 StopPoll() `void StopPoll ()`

10.116.3 Member Data Documentation

10.116.3.1 BOOST_BIT `const uint32_t BOOST_BIT = (1 << 0) [static]`

10.116.3.2 GND_SWITCH_BIT `const uint32_t GND_SWITCH_BIT = (1 << 1) [static]`

10.116.3.3 SYNC_BIT0 `const uint32_t SYNC_BIT0 = (1 << 2) [static]`


10.116.3.4 SYNC_BIT1 `const uint32_t SYNC_BIT1 = (1 << 3) [static]`

10.116.4 Event Documentation

10.116.4.1 PollStatusEvent `OnStgPollStatus^ PollStatusEvent`

10.117 CW2100DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CW2100DacqGroupChannelSelectionNet:

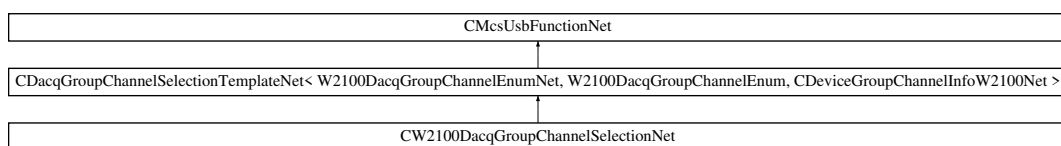

Public Member Functions

- `CW2100DacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)`

Additional Inherited Members


10.117.1 Constructor & Destructor Documentation

10.117.1.1 CW2100DacqGroupChannelSelectionNet() [CW2100DacqGroupChannelSelectionNet](#) (
[CMcsUsbNet](#)^ *mcsusb*)

10.118 CWarnerUssingDeviceNet Class Reference

[CWarnerUssingDeviceNet](#) is the class to control the Ussing device

Inheritance diagram for CWarnerUssingDeviceNet:

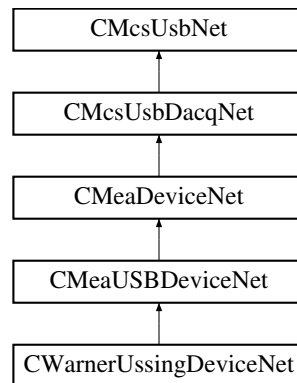

Public Member Functions

- [CWarnerUssingDeviceNet](#) ()
Initializes a new instance of the [CWarnerUssingDeviceNet](#) class.
- virtual [~CWarnerUssingDeviceNet](#) ()
- [ICWarnerUssingDeviceNet](#) ()

Properties

- [CWarnerUssingFunctionNet](#)^ [WarnerUssingFunction](#) [get]

Additional Inherited Members

10.118.1 Detailed Description

[CWarnerUssingDeviceNet](#) is the class to control the Ussing device

10.118.2 Constructor & Destructor Documentation

10.118.2.1 CWarnerUssingDeviceNet() CWarnerUssingDeviceNet ()

Initializes a new instance of the CWarnerUssingDeviceNet class.

10.118.2.2 ~CWarnerUssingDeviceNet() virtual ~CWarnerUssingDeviceNet () [virtual]

10.118.2.3 "!CWarnerUssingDeviceNet() !CWarnerUssingDeviceNet ()


10.118.3 Property Documentation

10.118.3.1 WarnerUssingFunction CWarnerUssingFunctionNet^ WarnerUssingFunction [get]

10.119 CWarnerUssingFunctionNet Class Reference

CWarnerUssingFunctionNet is the class to control the Ussing device

Inheritance diagram for CWarnerUssingFunctionNet:

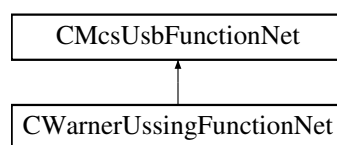

Public Member Functions

- [CWarnerUssingFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pWarner←
UssingFunctionPointerContainer)
Initializes a new instance of the [CWarnerUssingFunctionNet](#) class.
- [CWarnerUssingFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CWarnerUssingFunctionNet](#) ()
- [ICWarnerUssingFunctionNet](#) ()
- [int32_t](#) [GetChannelsCountOfChamber](#) ([int32_t](#) ChamberId)
gets number of channels in datastream from chamber amp with given index
- [int32_t](#) [GetNumberOfHardwareSlotsForChambers](#) ()
gets number of physical hardware slots for chambers amps
- [int32_t](#) [GetNumberOfAvailableChambers](#) ()
gets number of actually connected chamber amps
- [bool](#) [IsChamberAvailable](#) ([int32_t](#) ChamberId)
checks whether chamber amp is connected to slot
- void [SetPulse](#) ([int32_t](#) ChamberId, [UssingClampModeEnumNet](#) StgMode, [int32_t](#) NumberOfRepetitions,
array< [int](#) >[^] Amplitudes, array< [int](#) >[^] Durations)
defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes
- void [SetVoltageClampControllerParam_P](#) ([int32_t](#) ChamberId, [uint32_t](#) P)
sets P value of PID controller;
- void [SetVoltageClampControllerParam_I](#) ([int32_t](#) ChamberId, [uint32_t](#) I)
sets I value of PID controller;
- void [SetVoltageClampControllerParam_D](#) ([int32_t](#) ChamberId, [uint32_t](#) D)
sets D value of PID controller;
- [uint32_t](#) [GetVoltageClampControllerParam_P](#) ([int32_t](#) ChamberId)
gets P value of PID controller;
- [uint32_t](#) [GetVoltageClampControllerParam_I](#) ([int32_t](#) ChamberId)
gets I value of PID controller;
- [uint32_t](#) [GetVoltageClampControllerParam_D](#) ([int32_t](#) ChamberId)
gets D value of PID controller;
- void [SetClampMode](#) ([int32_t](#) ChamberId, [UssingClampModeEnumNet](#) NewClampMode)
sets clamp mode (voltage, current or open clamp)
- [UssingClampModeEnumNet](#) [GetClampMode](#) ([int32_t](#) ChamberId)
gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)
- [bool](#) [IsInternalCalibrationFinished](#) ([int32_t](#) ChamberId)
when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available
- [int32_t](#) [GetU1Offset](#) ([int32_t](#) ChamberId)
- [int32_t](#) [GetU2Offset](#) ([int32_t](#) ChamberId)
- [int32_t](#) [GetU1Reference](#) ([int32_t](#) ChamberId)
- [int32_t](#) [GetU2Reference](#) ([int32_t](#) ChamberId)
- [int32_t](#) [GetDacZero](#) ([int32_t](#) ChamberId)
- void [SetHighCurrentMode](#) ([int32_t](#) ChamberId)
- void [SetLowCurrentMode](#) ([int32_t](#) ChamberId)
- [bool](#) [IsHighCurrentMode](#) ([int32_t](#) ChamberId)
- [uint32_t](#) [GetLowCurrentRange](#) ([int32_t](#) ChamberId)
unit: nA
- [uint32_t](#) [GetHighCurrentRange](#) ([int32_t](#) ChamberId)
unit: nA
- [uint32_t](#) [GetDacPampsPerDigitLowCurrentRange](#) ([int32_t](#) ChamberId)
unit: pA/digit

- uint32_t [GetDacNampsPerDigitHighCurrentRange](#) (int32_t ChamberId)
unit: nA/digit
- uint32_t [GetUnitsPerDigit](#) (int32_t ChamberId, int32_t ChannelId)
gets amps/volts per digit for specified chamber and channel
- int32_t [GetUnitExponent](#) (int32_t ChamberId, int32_t ChannelId)
- UssingUnitEnumNet [GetUnitName](#) (int32_t ChamberId, int32_t ChannelId)
gets the channel's unit name
- String ^ [GetUnitDescription](#) (int32_t ChamberId, int32_t ChannelId)
- array< int > ^ [GetAvailableChambers](#) ()
returns array with (zero-based) ChamberIds of all available chambers
- int32_t [GetUptimeSeconds](#) (int32_t ChamberId)
- void [SetIdleModeOffset](#) (int32_t ChamberId, UssingClampModeEnumNet ClampMode, int32_t NewIdleOffset)
sets the offset (voltage or current) that will be applied when clamping is DISABLED
- int32_t [GetIdleModeOffset](#) (int32_t ChamberId, UssingClampModeEnumNet ClampMode)
gets the offset (voltage or current) that will be applied when clamping is DISABLED
- void [SetEnablePulse](#) (int32_t ChamberId, UssingClampModeEnumNet ClampMode, bool Enable)
enable pulse of given chamber and mode (voltage/current clamp) of this chamber
- bool [IsPulseEnabled](#) (int32_t ChamberId, UssingClampModeEnumNet ClampMode)
returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED
- void [SetLiquidResistance](#) (int32_t ChamberId, int32_t NewLiquidResistance_Ohm)
- int32_t [GetLiquidResistance](#) (int32_t ChamberId)
- int32_t [GetComplianceVoltage](#) (int32_t ChamberId)
returns compliance voltage threshold in digits; when Uc is above, current source is overloaded
- bool [CompensateElectrodeOffset](#) (int32_t ChamberId)
blocking call to compensate electrode offset of one chamber; returns true when successful
- bool [WaitForChamber](#) (int32_t ChamberId)
blocking call that waits for chamber boot-up calibration to complete
- bool [WaitForAllChambers](#) ()
blocking call that waits for ALL chambers' boot-up calibration to complete

Additional Inherited Members

10.119.1 Detailed Description

[CWarnerUssingFunctionNet](#) is the class to control the Ussing device

10.119.2 Constructor & Destructor Documentation

10.119.2.1 CWarnerUssingFunctionNet() [1/2] [CWarnerUssingFunctionNet](#) (
[CMcsUsbNet](#)^ mcsusb,
[CMcsUsbFunctionPointerContainer](#)^ pWarnerUssingFunctionPointerContainer)

Initializes a new instance of the [CWarnerUssingFunctionNet](#) class.

10.119.2.2 CWarnerUssingFunctionNet() [2/2] `CWarnerUssingFunctionNet (CMcsUsbNet^ mcsusb)`

10.119.2.3 ~CWarnerUssingFunctionNet() `virtual ~CWarnerUssingFunctionNet () [virtual]`

10.119.2.4 "!CWarnerUssingFunctionNet() `!CWarnerUssingFunctionNet ()`

10.119.3 Member Function Documentation

10.119.3.1 CompensateElectrodeOffset() `bool CompensateElectrodeOffset (int32_t ChamberId)`

blocking call to compensate electrode offset of one chamber; returns true when successful

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.2 GetAvailableChambers() `array<int> ^ GetAvailableChambers ()`

returns array with (zero-based) ChamberIds of all available chambers

10.119.3.3 GetChannelsCountOfChamber() `int32_t GetChannelsCountOfChamber (int32_t ChamberId)`

gets number of channels in datastream from chamber amp with given index

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

return value of zero means that amp is not placed

10.119.3.4 GetClampMode() `UssingClampModeEnumNet GetClampMode (`
`int32_t ChamberId)`

gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.5 GetComplianceVoltage() `int32_t GetComplianceVoltage (`
`int32_t ChamberId)`

returns compliance voltage threshold in digits; when U_c is above, current source is overloaded

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.6 GetDacNampsPerDigitHighCurrentRange() `uint32_t GetDacNampsPerDigitHighCurrentRange`
`(`
`int32_t ChamberId)`

unit: nA/digit

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.7 GetDacPampsPerDigitLowCurrentRange() `uint32_t GetDacPampsPerDigitLowCurrentRange (int32_t ChamberId)`

unit: pA/digit

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.8 GetDacZero() `int32_t GetDacZero (int32_t ChamberId)`

- diagnostic function only - ; gets real zero value of DAC in digits (0 -> neg. current; 32767 -> near zero; 65535 -> pos. current)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.9 GetHighCurrentRange() `uint32_t GetHighCurrentRange (int32_t ChamberId)`

unit: nA

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.10 GetIdleModeOffset() `int32_t GetIdleModeOffset (`
 `int32_t ChamberId,`
 `UssingClampModeEnumNet ClampMode)`

gets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation

Returns

unit: nA or mV

10.119.3.11 GetLiquidResistance() `int32_t GetLiquidResistance (`
 `int32_t ChamberId)`

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.12 GetLowCurrentRange() `uint32_t GetLowCurrentRange (`
 `int32_t ChamberId)`

unit: nA

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.13 GetNumberOfAvailableChambers() `int32_t GetNumberOfAvailableChambers ()`

gets number of actually connected chamber amps

Returns

10.119.3.14 GetNumberOfHardwareSlotsForChambers() `int32_t GetNumberOfHardwareSlotsForChambers ()`

gets number of physical hardware slots for chambers amps

Returns

10.119.3.15 GetU1Offset() `int32_t GetU1Offset (
int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.16 GetU1Reference() `int32_t GetU1Reference (
int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.17 GetU2Offset() `int32_t GetU2Offset (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.18 GetU2Reference() `int32_t GetU2Reference (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.19 GetUnitDescription() `String ^ GetUnitDescription (`
`int32_t ChamberId,`
`int32_t ChannelId)`

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

10.119.3.20 GetUnitExponent() `int32_t GetUnitExponent (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

example: return value -9 means that amps/volts per digit is in nano

10.119.3.21 GetUnitName() `UssingUnitEnumNet GetUnitName (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets the channel's unit name

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

10.119.3.22 GetUnitsPerDigit() `uint32_t GetUnitsPerDigit (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets amps/volts per digit for specified chamber and channel

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

10.119.3.23 GetUptimeSeconds() `int32_t GetUptimeSeconds (int32_t ChamberId)`

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.24 GetVoltageClampControllerParam_D() `uint32_t GetVoltageClampControllerParam_D (int32_t ChamberId)`

gets D value of PID controller;

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.25 GetVoltageClampControllerParam_I() `uint32_t GetVoltageClampControllerParam_I (int32_t ChamberId)`

gets I value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.26 GetVoltageClampControllerParam_P() `uint32_t GetVoltageClampControllerParam_P (`
`int32_t ChamberId)`

gets P value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.27 IsChamberAvailable() `bool IsChamberAvailable (`
`int32_t ChamberId)`

checks whether chamber amp is connected to slot

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

10.119.3.28 IsHighCurrentMode() `bool IsHighCurrentMode (`
`int32_t ChamberId)`

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.29 IsInternalCalibrationFinished() `bool IsInternalCalibrationFinished (
int32_t ChamberId)`

when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

10.119.3.30 IsPulseEnabled() `bool IsPulseEnabled (
int32_t ChamberId,
UssingClampModeEnumNet ClampMode)`

returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation

Returns

when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

10.119.3.31 SetClampMode() `void SetClampMode (
int32_t ChamberId,
UssingClampModeEnumNet NewClampMode)`

sets clamp mode (voltage, current or open clamp)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>NewClampMode</i>	

10.119.3.32 SetEnablePulse() void SetEnablePulse (
 int32_t *ChamberId*,
 UssingClampModeEnumNet *ClampMode*,
 bool *Enable*)

enable pulse of given chamber and mode (voltage/current clamp) of this chamber

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation
<i>Enable</i>	when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

10.119.3.33 SetHighCurrentMode() void SetHighCurrentMode (
 int32_t *ChamberId*)

Parameters

<i>ChamberId</i>	
------------------	--

10.119.3.34 SetIdleModeOffset() void SetIdleModeOffset (
 int32_t *ChamberId*,
 UssingClampModeEnumNet *ClampMode*,
 int32_t *NewIdleOffset*)

sets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation
<i>NewIdleOffset</i>	unit: nA or mV

10.119.3.35 SetLiquidResistance() void SetLiquidResistance (
 int32_t ChamberId,
 int32_t NewLiquidResistance_Ohm)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>NewLiquidResistance_Ohm</i>	

10.119.3.36 SetLowCurrentMode() void SetLowCurrentMode (
 int32_t ChamberId)

Parameters

<i>ChamberId</i>	
------------------	--

10.119.3.37 SetPulse() void SetPulse (
 int32_t ChamberId,
 UssingClampModeEnumNet StgMode,
 int32_t NumberOfRepetitions,
 array< int >^ Amplitudes,
 array< int >^ Durations)

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based); send pattern to connected amp
<i>StgMode</i>	voltage or current clamp stimulation
<i>NumberOfRepetitions</i>	number of repetitions for pulse pattern (-1 for infinite; n means pattern is applied n+1 times)
<i>Amplitudes</i>	amplitude; unit in voltage clamp: mV; unit in current clamp: nA
<i>Durations</i>	duration in 100us; CAUTION: first element is applied only one; auto-loop back to second element after last one

10.119.3.38 SetVoltageClampControllerParam_D() void SetVoltageClampControllerParam_D (
 int32_t ChamberId,
 uint32_t D)

sets D value of PID controller;

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>D</i>	useful range: 0..700

10.119.3.39 SetVoltageClampControllerParam_I() `void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)`

sets I value of PID controller;

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>I</i>	useful range: 80000..120000

10.119.3.40 SetVoltageClampControllerParam_P() `void SetVoltageClampControllerParam_P (int32_t ChamberId, uint32_t P)`

sets P value of PID controller;

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>P</i>	useful value: 130000

10.119.3.41 WaitForAllChambers() `bool WaitForAllChambers ()`

blocking call that waits for ALL chambers' boot-up calibration to complete

Returns

returns false when at least one chamber's calibration fails (e.g. timeout...)

10.119.3.42 WaitForChamber() `bool WaitForChamber (int32_t ChamberId)`

blocking call that waits for chamber boot-up calibration to complete

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---


Returns

returns false when calibration fails (e.g. timeout...)

10.120 CWarnerValveControllerDeviceNet Class Reference

[CWarnerValveControllerDeviceNet](#) is the class to access the Warner Valve Controller

Inheritance diagram for CWarnerValveControllerDeviceNet:

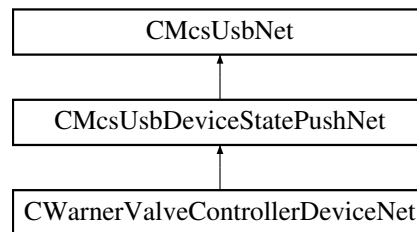

Public Member Functions

- delegate void [OnGetValveActive](#) (uint16_t valve, int valveActive)
- delegate void [OnGetValveManualState](#) (uint16_t valve, int32_t valveManualState)
- delegate void [OnGetValveManualGroup](#) (uint16_t valve, int32_t valveManualGroup)
- delegate void [OnGetValveMode](#) (uint16_t valve, WvcValveModeEnumNet ValveMode)
- delegate void [OnGetAnalogThresholdLow](#) (uint16_t valve, int32_t threshold)
- delegate void [OnGetAnalogThresholdHigh](#) (uint16_t valve, int32_t threshold)
- delegate void [OnGetDigitalPortDirection](#) (uint16_t port, PortDirectionEnumNet direction)
- delegate void [OnIsValveDigitalInInverted](#) (uint16_t valve, bool isInverted)
- delegate void [OnGetValveDigitalInPort](#) (uint16_t valve, uint32_t digitalInPort)
- delegate void [OnIsDigitalOutPortInverted](#) (uint16_t digitalOutPort, bool isInverted)
- delegate void [OnGetDigitalOutPortValve](#) (uint16_t digitalOutPort, uint32_t valve)
- delegate void [OnIsValveOpen](#) (uint16_t valve, bool valveOpen)
- delegate void [OnIsValveOpenInDigitalMode](#) (uint16_t valve, bool valveOpen)
- delegate void [OnIsValveOpenInAnalogMode](#) (uint16_t valve, bool valveOpen)
- delegate void [OnGetAnalogVoltage](#) (int32_t voltage)
- delegate void [OnTableEntryChanged](#) (uint16_t tableNumber)
- delegate void [OnGetTableNamebyIndex](#) (uint16_t tableNumber, String^ tableName)
- delegate void [OnGetActiveRunningTableNumber](#) (uint32_t tableNumber)
- delegate void [OnGetCurrentNumberOfValves](#) (int32_t numberOfValves)
- delegate void [OnGetValveBoardRevision](#) (uint32_t revision)
- delegate void [OnGetDisplayMode](#) (WvcDisplayModeEnumNet DisplayMode)
- [CWarnerValveControllerDeviceNet](#) ()
Initializes a new instance of the CWarnerValveControllerDeviceNet class.
- virtual [~CWarnerValveControllerDeviceNet](#) ()
- [!CWarnerValveControllerDeviceNet](#) ()
- int [GetValveActive](#) (uint16_t valve)

- Gets the valve active/inactive state*
- void [SetValveActive](#) (uint16_t valve, int valveActive)
- Sets the valve active/inactive state*
- uint32_t [GetValvesActiveMap](#) ()
- Gets the valves active/inactive states*
- void [SetValvesActiveMap](#) (uint32_t valvesActive)
- Sets the valve active/inactive state*
- int32_t [GetValveManualState](#) (uint16_t valve)
- Gets the valve manual on/off state*
- void [SetValveManualState](#) (uint16_t valve, int32_t valveManualState)
- Sets the valve manual on/off state*
- uint32_t [GetValvesManualStateMap](#) ()
- Gets the valves manual on/off states*
- void [SetValvesManualStateMap](#) (uint32_t valveaManualState)
- Sets the valve manual on/off state*
- int32_t [GetValveManualGroup](#) (uint16_t valve)
- Gets the valve manual group*
- void [SetValveManualGroup](#) (uint16_t valve, int32_t valveManualGroup)
- Sets the valve manual group*
- WvcValveModeEnumNet [GetValveMode](#) (uint16_t valve)
- Reads the valve mode*
- void [SetValveMode](#) (uint16_t valve, WvcValveModeEnumNet ValveMode)
- Writes the valve mode*
- int32_t [GetAnalogThresholdLow](#) (uint16_t valve)
- Gets the lower threshold for the analog in port per valve*
- void [SetAnalogThresholdLow](#) (uint16_t valve, int32_t threshold)
- Sets the lower threshold for the analog in port per valve*
- int32_t [GetAnalogThresholdHigh](#) (uint16_t valve)
- Gets the upper threshold for the analog in port per valve*
- void [SetAnalogThresholdHigh](#) (uint16_t valve, int32_t threshold)
- Sets the upper threshold for the analog in port per valve*
- PortDirectionEnumNet [GetDigitalPortDirection](#) (uint16_t port)
- Gets the direction of a digital port*
- void [SetDigitalPortDirection](#) (uint16_t port, PortDirectionEnumNet direction)
- Sets the direction of a digital port*
- bool [IsValveDigitalInInverted](#) (uint16_t valve)
- Is digital in inverted*
- void [SetValveDigitalInInvert](#) (uint16_t valve, bool isInverted)
- Invert digital in*
- uint32_t [GetValveDigitalInPort](#) (uint16_t valve)
- Gets the number of the digital in port which is mapped to a valve*
- void [SetValveDigitalInPort](#) (uint16_t valve, uint32_t digitalInPort)
- Map a digital in port to a valve*
- bool [IsDigitalOutPortInverted](#) (uint16_t digitalOutPort)
- Gets the number of the valve which is mapped to a digital out port*
- void [SetDigitalOutPortInvert](#) (uint16_t digitalOutPort, bool isInverted)
- Map a valve to a digital out port*
- uint32_t [GetDigitalOutPortValve](#) (uint16_t digitalOutPort)
- Gets the number of the valve which is mapped to a digital out port*
- void [SetDigitalOutPortValve](#) (uint16_t digitalOutPort, uint32_t valve)
- Map a valve to a digital out port*

- void [SetDefault](#) ()
Sets the settings of the valve controller to default
- bool [IsValveOpen](#) (uint16_t valve)
Is valve open
- bool [IsValveOpenInDigitalMode](#) (uint16_t valve)
True, if the valve would be open when the device is in digital mode
- bool [IsValveOpenInAnalogMode](#) (uint16_t valve)
True, if the valve would be open when the device is in analog mode
- int32_t [GetAnalogVoltage](#) ()
Reads the voltage on the analog in port
- void [GetValveTableEntry](#) (uint16_t valve, uint16_t index, [System::Runtime::InteropServices::Out]uint32_t% duration, [System::Runtime::InteropServices::Out]bool% state)
Read an entry from the valve protocol table
- void [SetValveTableEntry](#) (uint16_t valve, uint16_t index, uint32_t duration, bool state)
Write an entry to the valve protocol table
- void [ClearValveTable](#) (uint16_t valve)
Clear the valve protocol table
- void [LoadValveTable](#) ()
Load the current table from permanent memory
- void [StoreValveTable](#) ()
Store the current table in permanent memory
- String ^ [GetTableNamebyIndex](#) (uint16_t tableNumber)
Get the name of a protocol table
- String ^ [GetTableName](#) ()
Get the name of the current protocol table
- void [SetTableName](#) (String^ tableName)
Set the name of the current protocol table
- uint32_t [GetActiveRunningTableNumber](#) ()
Gets the number of the table that is active for running
- void [SetActiveRunningTableNumber](#) (uint32_t tableNumber)
Sets the number of the tanle that is active for running
- uint32_t [GetCurrentEditTableNumber](#) ()
Gets the number of the table that is current for editing
- void [SetCurrentEditTableNumber](#) (uint32_t tableNumber)
Sets the number of the table that is current for editing
- void [ClearTableName](#) ()
Clear the name of current protocol table
- void [SetTableStep](#) (uint16_t valve, int32_t steps)
Skips the table protocol for a valve by steps
- void [SetTableStepAll](#) (int32_t steps)
Skips the table protocol for all valves by steps
- int32_t [GetTotalNumberOfValves](#) ()
Get the total number of valves in the system
- int32_t [GetTotalNumberOfDigitalPorts](#) ()
Get the total number of digital ports in the system
- int32_t [GetTotalTableSize](#) ()
Get the total table size in the system
- int32_t [GetTotalNumberOfTables](#) ()
Get the total number of tables in the system
- int32_t [GetCurrentNumberOfValves](#) ()
Get the current number of valves connected to the system

- uint32_t [GetValveBoardRevision](#) ()
Gets the revision code of the valve board
- WvcDisplayModeEnumNet [GetDisplayMode](#) ()
Reads the display mode
- void [SetDisplayMode](#) (WvcDisplayModeEnumNet DisplayMode, int32_t lockTimeMs)
Writes the display mode
- String ^ [GetValveBoardRevisionString](#) ()
Gets the revision name of the valve board

Events

- [OnGetValveActive](#) ^ [GetValveActiveEvent](#) [add, remove, raise]
Event fires when the valve state for the valve number has changed
- [OnGetValveManualState](#) ^ [GetValveManualStateEvent](#) [add, remove, raise]
Event fires when the manual valve state for the valve number has changed
- [OnGetValveManualGroup](#) ^ [GetValveManualGroupEvent](#) [add, remove, raise]
Event fires when the manual valve group for the valve number has changed
- [OnGetValveMode](#) ^ [GetValveModeEvent](#) [add, remove, raise]
Event fires when the valve mode for the valve number has changed
- [OnGetAnalogThresholdLow](#) ^ [GetAnalogThresholdLowEvent](#) [add, remove, raise]
Event fires when the threshold in mV for the valve number has changed
- [OnGetAnalogThresholdHigh](#) ^ [GetAnalogThresholdHighEvent](#) [add, remove, raise]
Event fires when the threshold in mV for the valve number has changed
- [OnGetDigitalPortDirection](#) ^ [GetDigitalPortDirectionEvent](#) [add, remove, raise]
Event fires when the direction for the port number has changed
- [OnIsValveDigitalInInverted](#) ^ [IsValveDigitalInInvertedEvent](#) [add, remove, raise]
Event fires when is inverted for the valve number has changed
- [OnGetValveDigitalInPort](#) ^ [GetValveDigitalInPortEvent](#) [add, remove, raise]
Event fires when the digital in port for the valve number has changed
- [OnIsDigitalOutPortInverted](#) ^ [IsDigitalOutPortInvertedEvent](#) [add, remove, raise]
Event fires when is inverted for the digital out port has changed
- [OnGetDigitalOutPortValve](#) ^ [GetDigitalOutPortValveEvent](#) [add, remove, raise]
Event fires when the valve number for the digital out port has changed
- [OnIsValveOpen](#) ^ [IsValveOpenEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnIsValveOpenInDigitalMode](#) ^ [IsValveOpenInDigitalModeEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnIsValveOpenInAnalogMode](#) ^ [IsValveOpenInAnalogModeEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnGetAnalogVoltage](#) ^ [GetAnalogVoltageEvent](#) [add, remove, raise]
Event fires when the voltage in mV has changed
- [OnTableEntryChanged](#) ^ [TableEntryChangedEvent](#) [add, remove, raise]
Event fires when an entry of a table changed
- [OnGetTableNamebyIndex](#) ^ [GetTableNamebyIndexEvent](#) [add, remove, raise]
Event fires when the name of the table for the table number has changed
- [OnGetActiveRunningTableNumber](#) ^ [GetActiveRunningTableNumberEvent](#) [add, remove, raise]
Event fires when the table number has changed
- [OnGetCurrentNumberOfValves](#) ^ [GetCurrentNumberOfValvesEvent](#) [add, remove, raise]
Event fires when the number of valves has changed
- [OnGetValveBoardRevision](#) ^ [GetValveBoardRevisionEvent](#) [add, remove, raise]
Event fires when the revision code has changed
- [OnGetDisplayMode](#) ^ [GetDisplayModeEvent](#) [add, remove, raise]
Event fires when the display mode has changed

Additional Inherited Members

10.120.1 Detailed Description

[CWarnerValveControllerDeviceNet](#) is the class to access the Warner Valve Controller

10.120.2 Constructor & Destructor Documentation

10.120.2.1 CWarnerValveControllerDeviceNet() [CWarnerValveControllerDeviceNet](#) ()

Initializes a new instance of the [CWarnerValveControllerDeviceNet](#) class.

10.120.2.2 ~CWarnerValveControllerDeviceNet() [virtual](#) [~CWarnerValveControllerDeviceNet](#) () [virtual]

10.120.2.3 ~!CWarnerValveControllerDeviceNet() [!CWarnerValveControllerDeviceNet](#) ()

10.120.3 Member Function Documentation

10.120.3.1 ClearTableName() [void](#) [ClearTableName](#) ()

Clear the name of current protocol table

10.120.3.2 ClearValveTable() [void](#) [ClearValveTable](#) ([uint16_t](#) *valve*)

Clear the valve protocol table

Parameters

<i>valve</i>	The valve number
--------------	------------------

10.120.3.3 GetActiveRunningTableNumber() [uint32_t](#) [GetActiveRunningTableNumber](#) ()

Gets the number of the table that is active for running

Returns

The table number

10.120.3.4 GetAnalogThresholdHigh() `int32_t GetAnalogThresholdHigh (uint16_t valve)`

Gets the upper threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The threshold in mV

10.120.3.5 GetAnalogThresholdLow() `int32_t GetAnalogThresholdLow (uint16_t valve)`

Gets the lower threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The threshold in mV

10.120.3.6 GetAnalogVoltage() `int32_t GetAnalogVoltage ()`

Reads the voltage on the analog in port

Returns

The voltage in mV

10.120.3.7 GetCurrentEditTableNumber() `uint32_t GetCurrentEditTableNumber ()`

Gets the number of the table that is current for editing

Returns

The table number

10.120.3.8 GetCurrentNumberOfValves() `int32_t GetCurrentNumberOfValves ()`

Get the current number of valves connected to the system

Returns

The number of valves

10.120.3.9 GetDigitalOutPortValve() `uint32_t GetDigitalOutPortValve (
uint16_t digitalOutPort)`

Gets the number of the valve which is mapped to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
-----------------------	----------------------

Returns

The valve number

10.120.3.10 GetDigitalPortDirection() `PortDirectionEnumNet GetDigitalPortDirection (
uint16_t port)`

Gets the direction of a digital port

Parameters

<i>port</i>	The port number
-------------	-----------------

Returns

the direction

10.120.3.11 GetDisplayMode() `WvcDisplayModeEnumNet GetDisplayMode ()`

Reads the display mode

Returns

The display mode

10.120.3.12 GetTableName() `String ^ GetTableName ()`

Get the name of the current protocol table

Returns

The name of the table

10.120.3.13 GetTableNamebyIndex() `String ^ GetTableNamebyIndex (
uint16_t tableNumber)`

Get the name of a protocol table

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

Returns

The name of the table

10.120.3.14 GetTotalNumberOfDigitalPorts() `int32_t GetTotalNumberOfDigitalPorts ()`

Get the total number of digital ports in the system

Returns

The number of digital ports

10.120.3.15 GetTotalNumberOfTables() `int32_t GetTotalNumberOfTables ()`

Get the total number of tables in the system

Returns

The number of tables

10.120.3.16 GetTotalNumberOfValves() `int32_t GetTotalNumberOfValves ()`

Get the total number of valves in the system

Returns

The number of valves

10.120.3.17 GetTotalTableSize() `int32_t GetTotalTableSize ()`

Get the total table size in the system

Returns

The table size

10.120.3.18 GetValveActive() `int GetValveActive (uint16_t valve)`

Gets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve state

10.120.3.19 GetValveBoardRevision() `uint32_t GetValveBoardRevision ()`

Gets the revision code of the valve board

Returns

The revision code

10.120.3.20 GetValveBoardRevisionString() `String ^ GetValveBoardRevisionString ()`

Gets the revision name of the valve board

Returns

The revision name

10.120.3.21 GetValveDigitalInPort() `uint32_t GetValveDigitalInPort (uint16_t valve)`

Gets the number of the digital in port which is mapped to a valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The digital in port

10.120.3.22 GetValveManualGroup() `int32_t GetValveManualGroup (uint16_t valve)`

Gets the valve manual group

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The manual valve group

10.120.3.23 GetValveManualState() `int32_t GetValveManualState (uint16_t valve)`

Gets the valve manual on/off state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The manual valve state

10.120.3.24 GetValveMode() `WvcValveModeEnumNet GetValveMode (uint16_t valve)`

Reads the valve mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve mode

10.120.3.25 GetValvesActiveMap() `uint32_t GetValvesActiveMap ()`

Gets the valves active/inactive states

Returns

The valves states

10.120.3.26 GetValvesManualStateMap() `uint32_t GetValvesManualStateMap ()`

Gets the valves manual on/off states

Returns

The manual valves states

10.120.3.27 GetValveTableEntry() `void GetValveTableEntry (
 uint16_t valve,
 uint16_t index,
 [System::Runtime::InteropServices::Out] uint32_t% duration,
 [System::Runtime::InteropServices::Out] bool% state)`

Read an entry from the valve protocol table

Parameters

<i>valve</i>	The valve number
<i>index</i>	The index in the table
<i>duration</i>	the duration in ms
<i>state</i>	the state

10.120.3.28 IsDigitalOutPortInverted() `bool IsDigitalOutPortInverted (
 uint16_t digitalOutPort)`

Gets the number of the valve which is mapped to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
-----------------------	----------------------

Returns

is inverted

10.120.3.29 IsValveDigitalInInverted() `bool IsValveDigitalInInverted (uint16_t valve)`

Is digital in inverted

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is inverted

10.120.3.30 IsValveOpen() `bool IsValveOpen (uint16_t valve)`

Is valve open

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

10.120.3.31 IsValveOpenInAnalogMode() `bool IsValveOpenInAnalogMode (uint16_t valve)`

True, if the valve would be open when the device is in analog mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

10.120.3.32 IsValveOpenInDigitalMode() `bool IsValveOpenInDigitalMode (uint16_t valve)`

True, if the valve would be open when the device is in digital mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

10.120.3.33 LoadValveTable() `void LoadValveTable ()`

Load the current table from permanent memory

10.120.3.34 OnGetActiveRunningTableNumber() `delegate void OnGetActiveRunningTableNumber (uint32_t tableNumber)`

10.120.3.35 OnGetAnalogThresholdHigh() `delegate void OnGetAnalogThresholdHigh (uint16_t valve, int32_t threshold)`

10.120.3.36 OnGetAnalogThresholdLow() `delegate void OnGetAnalogThresholdLow (uint16_t valve, int32_t threshold)`

10.120.3.37 OnGetAnalogVoltage() delegate void OnGetAnalogVoltage (
int32_t voltage)

10.120.3.38 OnGetCurrentNumberOfValves() delegate void OnGetCurrentNumberOfValves (
int32_t numberOfValves)

10.120.3.39 OnGetDigitalOutPortValve() delegate void OnGetDigitalOutPortValve (
uint16_t digitalOutPort,
uint32_t valve)

10.120.3.40 OnGetDigitalPortDirection() delegate void OnGetDigitalPortDirection (
uint16_t port,
PortDirectionEnumNet direction)

10.120.3.41 OnGetDisplayMode() delegate void OnGetDisplayMode (
WvcDisplayModeEnumNet DisplayMode)

10.120.3.42 OnGetTableNamebyIndex() delegate void OnGetTableNamebyIndex (
uint16_t tableNumber,
String^ tableName)

10.120.3.43 OnGetValveActive() delegate void OnGetValveActive (
uint16_t valve,
int valveActive)

10.120.3.44 OnGetValveBoardRevision() delegate void OnGetValveBoardRevision (
uint32_t revision)

10.120.3.45 OnGetValveDigitalInPort() delegate void OnGetValveDigitalInPort (
uint16_t valve,
uint32_t digitalInPort)

10.120.3.46 OnGetValveManualGroup() delegate void OnGetValveManualGroup (
 uint16_t valve,
 int32_t valveManualGroup)

10.120.3.47 OnGetValveManualState() delegate void OnGetValveManualState (
 uint16_t valve,
 int32_t valveManualState)

10.120.3.48 OnGetValveMode() delegate void OnGetValveMode (
 uint16_t valve,
 WvcValveModeEnumNet ValveMode)

10.120.3.49 OnIsDigitalOutPortInverted() delegate void OnIsDigitalOutPortInverted (
 uint16_t digitalOutPort,
 bool isInverted)

10.120.3.50 OnIsValveDigitalInInverted() delegate void OnIsValveDigitalInInverted (
 uint16_t valve,
 bool isInverted)

10.120.3.51 OnIsValveOpen() delegate void OnIsValveOpen (
 uint16_t valve,
 bool valveOpen)

10.120.3.52 OnIsValveOpenInAnalogMode() delegate void OnIsValveOpenInAnalogMode (
 uint16_t valve,
 bool valveOpen)

10.120.3.53 OnIsValveOpenInDigitalMode() delegate void OnIsValveOpenInDigitalMode (
 uint16_t valve,
 bool valveOpen)

10.120.3.54 OnTableEntryChanged() delegate void OnTableEntryChanged (
 uint16_t tableNumber)

10.120.3.55 SetActiveRunningTableNumber() void SetActiveRunningTableNumber (
 uint32_t tableNumber)

Sets the number of the tanle that is active for running

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

10.120.3.56 SetAnalogThresholdHigh() `void SetAnalogThresholdHigh (`
 `uint16_t valve,`
 `int32_t threshold)`

Sets the upper threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
<i>threshold</i>	The threshold in mV

10.120.3.57 SetAnalogThresholdLow() `void SetAnalogThresholdLow (`
 `uint16_t valve,`
 `int32_t threshold)`

Sets the lower threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
<i>threshold</i>	The threshold in mV

10.120.3.58 SetCurrentEditTableNumber() `void SetCurrentEditTableNumber (`
 `uint32_t tableNumber)`

Sets the number of the table that is current for editing

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

10.120.3.59 SetDefault() `void SetDefault ()`

Sets the settings of the valve controller to default

10.120.3.60 SetDigitalOutPortInvert() void SetDigitalOutPortInvert (
 uint16_t *digitalOutPort*,
 bool *isInverted*)

Map a valve to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
<i>isInverted</i>	True if digital out is to be inverted

10.120.3.61 SetDigitalOutPortValve() void SetDigitalOutPortValve (
 uint16_t *digitalOutPort*,
 uint32_t *valve*)

Map a valve to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
<i>valve</i>	The valve number

10.120.3.62 SetDigitalPortDirection() void SetDigitalPortDirection (
 uint16_t *port*,
 PortDirectionEnumNet *direction*)

Sets the direction of a digital port

Parameters

<i>port</i>	The port number
<i>direction</i>	the direction

10.120.3.63 SetDisplayMode() void SetDisplayMode (
 WvcDisplayModeEnumNet *DisplayMode*,
 int32_t *lockTimeMs*)

Writes the display mode

Parameters

<i>DisplayMode</i>	The display mode
<i>lockTimeMs</i>	Locks the display for ms

10.120.3.64 SetTableName() void SetTableName (
String^ *tableName*)

Set the name of the current protocol table

Parameters

<i>tableName</i>	The name of the table
------------------	-----------------------

10.120.3.65 SetTableStep() void SetTableStep (
uint16_t *valve*,
int32_t *steps*)

Skips the table protocol for a valve by steps

Parameters

<i>valve</i>	The valve number
<i>steps</i>	Number of steps

10.120.3.66 SetTableStepAll() void SetTableStepAll (
int32_t *steps*)

Skips the table protocol for all valves by steps

Parameters

<i>steps</i>	Number of steps
--------------	-----------------

10.120.3.67 SetValveActive() void SetValveActive (
uint16_t *valve*,
int *valveActive*)

Sets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
<i>valveActive</i>	The valve state

10.120.3.68 SetValveDigitalInInvert() `void SetValveDigitalInInvert (`
 `uint16_t valve,`
 `bool isInverted)`

Invert digital in

Parameters

<i>valve</i>	The valve number
<i>isInverted</i>	True if digital in is to be inverted

10.120.3.69 SetValveDigitalInPort() `void SetValveDigitalInPort (`
 `uint16_t valve,`
 `uint32_t digitalInPort)`

Map a digital in port to a valve

Parameters

<i>valve</i>	The valve number
<i>digitalInPort</i>	The digital in port

10.120.3.70 SetValveManualGroup() `void SetValveManualGroup (`
 `uint16_t valve,`
 `int32_t valveManualGroup)`

Sets the valve manual group

Parameters

<i>valve</i>	The valve number
<i>valveManualGroup</i>	The manual valve group

10.120.3.71 SetValveManualState() `void SetValveManualState (`
 `uint16_t valve,`
 `int32_t valveManualState)`

Sets the valve manual on/off state

Parameters

<i>valve</i>	The valve number
<i>valveManualState</i>	The manual valve state

10.120.3.72 SetValveMode() void SetValveMode (
 uint16_t *valve*,
 WvcValveModeEnumNet *ValveMode*)

Writes the valve mode

Parameters

<i>valve</i>	The valve number
<i>ValveMode</i>	The valve mode

10.120.3.73 SetValvesActiveMap() void SetValvesActiveMap (
 uint32_t *valvesActive*)

Sets the valve active/inactive state

Parameters

<i>valvesActive</i>	The valves states
---------------------	-------------------

10.120.3.74 SetValvesManualStateMap() void SetValvesManualStateMap (
 uint32_t *valveaManualState*)

Sets the valve manual on/off state

Parameters

<i>valveaManualState</i>	The manual valves states
--------------------------	--------------------------

10.120.3.75 SetValveTableEntry() void SetValveTableEntry (
 uint16_t *valve*,
 uint16_t *index*,
 uint32_t *duration*,
 bool *state*)

Write an entry to the valve protocol table

Parameters

<i>valve</i>	The valve number
<i>index</i>	The index in the table
<i>duration</i>	the duration in ms
<i>state</i>	the state

10.120.3.76 StoreValveTable() `void StoreValveTable ()`

Store the current table in permanent memory

10.120.4 Event Documentation**10.120.4.1 GetActiveRunningTableNumberEvent** `OnGetActiveRunningTableNumber^ GetActiveRunning↔
TableNumberEvent [add], [remove], [raise]`

Event fires when the table number has changed

10.120.4.2 GetAnalogThresholdHighEvent `OnGetAnalogThresholdHigh^ GetAnalogThresholdHighEvent
[add], [remove], [raise]`

Event fires when the threshold in mV for the valve number has changed

10.120.4.3 GetAnalogThresholdLowEvent `OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent
[add], [remove], [raise]`

Event fires when the threshold in mV for the valve number has changed

10.120.4.4 GetAnalogVoltageEvent `OnGetAnalogVoltage^ GetAnalogVoltageEvent [add], [remove],
[raise]`

Event fires when the voltage in mV has changed

10.120.4.5 GetCurrentNumberOfValvesEvent [OnGetCurrentNumberOfValves](#)[^] [GetCurrentNumberOfValvesEvent](#) [add], [remove], [raise]

Event fires when the number of valves has changed

10.120.4.6 GetDigitalOutPortValveEvent [OnGetDigitalOutPortValve](#)[^] [GetDigitalOutPortValveEvent](#) [add], [remove], [raise]

Event fires when the valve number for the digital out port has changed

10.120.4.7 GetDigitalPortDirectionEvent [OnGetDigitalPortDirection](#)[^] [GetDigitalPortDirectionEvent](#) [add], [remove], [raise]

Event fires when the direction for the port number has changed

10.120.4.8 GetDisplayModeEvent [OnGetDisplayMode](#)[^] [GetDisplayModeEvent](#) [add], [remove], [raise]

Event fires when the display mode has changed

10.120.4.9 GetTableNamebyIndexEvent [OnGetTableNamebyIndex](#)[^] [GetTableNamebyIndexEvent](#) [add], [remove], [raise]

Event fires when the name of the table for the table number has changed

10.120.4.10 GetValveActiveEvent [OnGetValveActive](#)[^] [GetValveActiveEvent](#) [add], [remove], [raise]

Event fires when the valve state for the valve number has changed

10.120.4.11 GetValveBoardRevisionEvent [OnGetValveBoardRevision](#)[^] [GetValveBoardRevisionEvent](#) [add], [remove], [raise]

Event fires when the revision code has changed

10.120.4.12 GetValveDigitalInPortEvent [OnGetValveDigitalInPort](#)[^] GetValveDigitalInPortEvent [add], [remove], [raise]

Event fires when the digital in port for the valve number has changed

10.120.4.13 GetValveManualGroupEvent [OnGetValveManualGroup](#)[^] GetValveManualGroupEvent [add], [remove], [raise]

Event fires when the manual valve group for the valve number has changed

10.120.4.14 GetValveManualStateEvent [OnGetValveManualState](#)[^] GetValveManualStateEvent [add], [remove], [raise]

Event fires when the manual valve state for the valve number has changed

10.120.4.15 GetValveModeEvent [OnGetValveMode](#)[^] GetValveModeEvent [add], [remove], [raise]

Event fires when the valve mode for the valve number has changed

10.120.4.16 IsDigitalOutPortInvertedEvent [OnIsDigitalOutPortInverted](#)[^] IsDigitalOutPortInverted↔ Event [add], [remove], [raise]

Event fires when is inverted for the digital out port has changed

10.120.4.17 IsValveDigitalInInvertedEvent [OnIsValveDigitalInInverted](#)[^] IsValveDigitalInInverted↔ Event [add], [remove], [raise]

Event fires when is inverted for the valve number has changed

10.120.4.18 IsValveOpenEvent [OnIsValveOpen](#)[^] IsValveOpenEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed

10.120.4.19 IsValveOpenInAnalogModeEvent [OnIsValveOpenInAnalogMode](#)[^] [IsValveOpenInAnalogMode](#)↔
Event [add], [remove], [raise]

Event fires when is open for the valve number has changed

10.120.4.20 IsValveOpenInDigitalModeEvent [OnIsValveOpenInDigitalMode](#)[^] [IsValveOpenInDigital](#)↔
ModeEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed


10.120.4.21 TableEntryChangedEvent [OnTableEntryChanged](#)[^] [TableEntryChangedEvent](#) [add], [remove], [raise]

Event fires when an entry of a table changed

10.121 CWarnerValveControllerDeviceTesterFunctionNet Class Reference

[CWarnerValveControllerDeviceTesterFunctionNet](#) is the class to access the functions for the Warner Valve Controller Device Tester

Inheritance diagram for CWarnerValveControllerDeviceTesterFunctionNet:

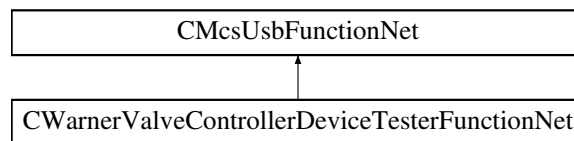

Public Member Functions

- [CWarnerValveControllerDeviceTesterFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pWarnerValveControllerDeviceTesterFunctionPointerContainer)
Initializes a new instance of the [CWarnerValveControllerDeviceTesterFunctionNet](#) class.
- [CWarnerValveControllerDeviceTesterFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CWarnerValveControllerDeviceTesterFunctionNet](#) ()
- [!CWarnerValveControllerDeviceTesterFunctionNet](#) ()
- void [SetADC](#) (uint32_t onoff)
Sets the ADC port of the tester
- uint32_t [GetSync](#) ()
Gets the output from the sync port
- void [SetTrigger](#) (uint32_t trigger)
Sets the input to the trigger port
- void [SetTriggerSyncDirection](#) (uint32_t direction)
Sets the direction of the trigger/sync test port
- uint32_t [GetIO](#) ()
Gets the output from the io ports
- void [SetIO](#) (uint32_t io)
Sets the input to the io ports
- void [SetIODirection](#) (uint32_t direction)
Sets the direction of the IO test ports

Additional Inherited Members

10.121.1 Detailed Description

[CWarnerValveControllerDeviceTesterFunctionNet](#) is the class to access the functions for the Warner Valve Controller Device Tester

10.121.2 Constructor & Destructor Documentation

10.121.2.1 CWarnerValveControllerDeviceTesterFunctionNet() [1/2] [CWarnerValveControllerDeviceTesterFunctionNet](#)
(
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pWarnerValveControllerDeviceTesterFunction*↔
 PointerContainer)

Initializes a new instance of the [CWarnerValveControllerDeviceTesterFunctionNet](#) class.

10.121.2.2 CWarnerValveControllerDeviceTesterFunctionNet() [2/2] [CWarnerValveControllerDeviceTesterFunctionNet](#)
(
 [CMcsUsbNet](#)[^] *mcsusb*)

10.121.2.3 ~CWarnerValveControllerDeviceTesterFunctionNet() virtual [~CWarnerValveControllerDeviceTesterFunctionNet](#)
() [virtual]

10.121.2.4 "!CWarnerValveControllerDeviceTesterFunctionNet() [!CWarnerValveControllerDeviceTesterFunctionNet](#)
()

10.121.3 Member Function Documentation

10.121.3.1 GetIO() uint32_t GetIO ()

Gets the output from the io ports

Returns

The manual valves states

10.121.3.2 GetSync() `uint32_t GetSync ()`

Gets the output from the sync port

Returns

The sync state

10.121.3.3 SetADC() `void SetADC (
uint32_t onoff)`

Sets the ADC port of the tester

Parameters

<i>onoff</i>	The port state
--------------	----------------

10.121.3.4 SetIO() `void SetIO (
uint32_t io)`

Sets the input to the io ports

Parameters

<i>io</i>	The manual valves states
-----------	--------------------------

10.121.3.5 SetIODirection() `void SetIODirection (
int32_t direction)`

Sets the direction of the IO test ports

Parameters

<i>direction</i>	The 16bit direction map: 1=IN 0=OUT
------------------	-------------------------------------

10.121.3.6 SetTrigger() `void SetTrigger (
uint32_t trigger)`

Sets the input to the trigger port

Parameters

<i>trigger</i>	The trigger state
----------------	-------------------

10.121.3.7 SetTriggerSyncDirection() `void SetTriggerSyncDirection (uint32_t direction)`


Sets the direction of the trigger/sync test port

Parameters

<i>direction</i>	The direction: 1=IN 0=OUT
------------------	---------------------------

10.122 CWClassicFunctionNet Class Reference

Inheritance diagram for CWClassicFunctionNet:

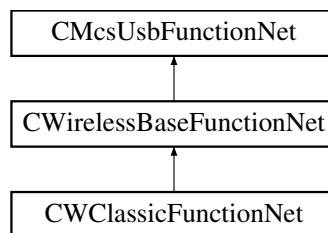

Public Member Functions

- [CWClassicFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] wClassicFunctionPointerContainer)
- [CWClassicFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- `uint32_t` [ResetChannelmap](#) (unsigned int virtualDevice)
- `uint32_t` [SetChannelmap](#) (unsigned char position, unsigned char channel, unsigned int Device)
- `void` [SetHWSelectedChannels](#) (array< bool >[^] channels, unsigned int Device)
- `void` [SetRFLostBehaviour](#) (uint8_t stoponfailure, unsigned int Device)
- `void` [SetHeadstageOnOff](#) (uint16_t onoff)
- `USHORT` [GetHeadstageOnOff](#) ()
- `void` [SetRFFrequencyHeadstage](#) (uint8_t receiver_nb, unsigned short frequency)
- `unsigned short` [GetRFFrequencyHeadstage](#) (uint8_t receiver_nb)
- `void` [SetRFFrequencyReceiver](#) (uint8_t receiver_nb, uint8_t configuration, unsigned short frequency)
- `void` [SetRFFrequencyReceiverEeprom](#) (uint8_t receiver_nb, uint8_t configuration, unsigned short frequency)
- `unsigned short` [GetRFFrequencyReceiver](#) (uint8_t receiver_nb, uint8_t configuration)
- `void` [SetSerialNumberHeadstage](#) (unsigned short number)
- `unsigned short` [GetSerialNumberHeadstage](#) ()
- `void` [SetSelectedHeadstage](#) (uint8_t number)
- `uint8_t` [GetSelectedHeadstage](#) ()
- `void` [ScanForHeadstages](#) ()
- `uint8_t` [GetScanHeadstagesResult](#) (int max_wait_for_ms)

- void [SetFilterParametersHeadstage](#) (unsigned short index, array< int > ^ buffer)
- array< int > ^ [GetFilterParametersHeadstage](#) (unsigned short index)
- bool [GetHasRedLedHeadstage](#) ()
- void [SetHasChecksum](#) (unsigned int has, unsigned int Device)
- unsigned int [GetHasChecksum](#) (unsigned int Device)
- void [SetResetFilter](#) (unsigned int reset, unsigned int Device)
- unsigned int [GetResetFilter](#) (unsigned int Device)
- void [SetWPAType](#) (unsigned short type, unsigned int Device)
- unsigned short [GetWPAType](#) (unsigned int Device)
- void [SetWPADebugMode](#) (unsigned int mode, unsigned int Device)
- unsigned int [GetWPADebugMode](#) (unsigned int Device)
- void [SetRFPower](#) (unsigned short power)
- unsigned short [GetRFPower](#) ()
- unsigned int [GetRFConnectionStatus](#) ()

Additional Inherited Members

10.122.1 Constructor & Destructor Documentation

10.122.1.1 CWClassicFunctionNet() [1/2] [CWClassicFunctionNet](#) (
[CMcsUsbNet](#) ^ *mcsusb*,
[CMcsUsbFunctionPointerContainer](#) ^ *wClassicFuntionPointerContainer*)

10.122.1.2 CWClassicFunctionNet() [2/2] [CWClassicFunctionNet](#) (
[CMcsUsbNet](#) ^ *mcsusb*)

10.122.2 Member Function Documentation

10.122.2.1 GetFilterParametersHeadstage() array<int> ^ [GetFilterParametersHeadstage](#) (
 unsigned short *index*)

10.122.2.2 GetHasChecksum() unsigned int [GetHasChecksum](#) (
 unsigned int *Device*)

10.122.2.3 GetHasRedLedHeadstage() bool [GetHasRedLedHeadstage](#) ()

10.122.2.4 GetHeadstageOnOff() `USHORT GetHeadstageOnOff ()`

10.122.2.5 GetResetFilter() `unsigned int GetResetFilter (`
`unsigned int Device)`

10.122.2.6 GetRFConnectionStatus() `unsigned int GetRFConnectionStatus ()`

10.122.2.7 GetRFFrequencyHeadstage() `unsigned short GetRFFrequencyHeadstage (`
`uint8_t receiver_nb)`

10.122.2.8 GetRFFrequencyReceiver() `unsigned short GetRFFrequencyReceiver (`
`uint8_t receiver_nb,`
`uint8_t configuration)`

10.122.2.9 GetRFPower() `unsigned short GetRFPower ()`

10.122.2.10 GetScanHeadstagesResult() `uint8_t GetScanHeadstagesResult (`
`int max_wait_for_ms)`

10.122.2.11 GetSelectedHeadstage() `uint8_t GetSelectedHeadstage ()`

10.122.2.12 GetSerialNumberHeadstage() `unsigned short GetSerialNumberHeadstage ()`

10.122.2.13 GetWPADebugMode() `unsigned int GetWPADebugMode (`
`unsigned int Device)`

10.122.2.14 GetWPAType() unsigned short GetWPAType (
 unsigned int *Device*)

10.122.2.15 ResetChannelmap() uint32_t ResetChannelmap (
 unsigned int *virtualDevice*)

10.122.2.16 ScanForHeadstages() void ScanForHeadstages ()

10.122.2.17 SetChannelmap() uint32_t SetChannelmap (
 unsigned char *position*,
 unsigned char *channel*,
 unsigned int *Device*)

10.122.2.18 SetFilterParametersHeadstage() void SetFilterParametersHeadstage (
 unsigned short *index*,
 array< int >^ *buffer*)

10.122.2.19 SetHasChecksum() void SetHasChecksum (
 unsigned int *has*,
 unsigned int *Device*)

10.122.2.20 SetHeadstageOnOff() void SetHeadstageOnOff (
 uint16_t *onoff*)


10.122.2.21 SetHWSelectedChannels() void SetHWSelectedChannels (
 array< bool >^ *channels*,
 unsigned int *Device*)

10.122.2.22 SetResetFilter() void SetResetFilter (
 unsigned int *reset*,
 unsigned int *Device*)

- 10.122.2.23 SetRFFrequencyHeadstage()** void SetRFFrequencyHeadstage (
 uint8_t *receiver_nb*,
 unsigned short *frequency*)
- 10.122.2.24 SetRFFrequencyReceiver()** void SetRFFrequencyReceiver (
 uint8_t *receiver_nb*,
 uint8_t *configuration*,
 unsigned short *frequency*)
- 10.122.2.25 SetRFFrequencyReceiverEeprom()** void SetRFFrequencyReceiverEeprom (
 uint8_t *receiver_nb*,
 uint8_t *configuration*,
 unsigned short *frequency*)
- 10.122.2.26 SetRFLostBehaviour()** void SetRFLostBehaviour (
 uint8_t *stoponfailure*,
 unsigned int *Device*)
- 10.122.2.27 SetRFPower()** void SetRFPower (
 unsigned short *power*)
- 10.122.2.28 SetSelectedHeadstage()** void SetSelectedHeadstage (
 uint8_t *number*)
- 10.122.2.29 SetSerialNumberHeadstage()** void SetSerialNumberHeadstage (
 unsigned short *number*)
- 10.122.2.30 SetWPADebugMode()** void SetWPADebugMode (
 unsigned int *mode*,
 unsigned int *Device*)
- 10.122.2.31 SetWPAType()** void SetWPAType (
 unsigned short *type*,
 unsigned int *Device*)

10.123 CWirelessBaseFunctionNet Class Reference

Inheritance diagram for CWirelessBaseFunctionNet:

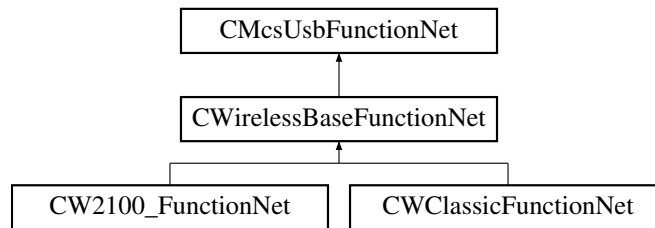

Public Member Functions

- [CWirelessBaseFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] mcsusbfunction)

Static Public Member Functions

- static String [^] [CreateWirelessHeadstageSerialNumberString](#) (unsigned short ID)

Additional Inherited Members

10.123.1 Constructor & Destructor Documentation

10.123.1.1 CWirelessBaseFunctionNet() [CWirelessBaseFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] mcsusbfunction)

10.123.2 Member Function Documentation

10.123.2.1 CreateWirelessHeadstageSerialNumberString() static String [^] [CreateWirelessHeadstageSerialNumberString](#) (
 unsigned short ID) [static]

10.124 DeviceIdNet Struct Reference

Device Id.

Public Member Functions

- [DeviceIdNet](#) ()
- [DeviceIdNet](#) (VendorIdEnumNet vendor, ProductIdEnumNet product, int bcd, McsBusTypeEnumNet bustype)
- [DeviceIdNet](#) ([DeviceIdNet](#)% deviceId)
- [DeviceIdNet](#) operator= ([DeviceIdNet](#)% deviceId)

Public Attributes

- VendorIdEnumNet [IdVendor](#)
- ProductIdEnumNet [IdProduct](#)
- int [BcdDevice](#)
- McsBusTypeEnumNet [BusType](#)

10.124.1 Detailed Description

Device Id.

10.124.2 Constructor & Destructor Documentation

10.124.2.1 DeviceIdNet() [1/3] [DeviceIdNet](#) ()

10.124.2.2 DeviceIdNet() [2/3] [DeviceIdNet](#) (VendorIdEnumNet *vendor*, ProductIdEnumNet *product*, int *bcd*, McsBusTypeEnumNet *bustype*)

10.124.2.3 DeviceIdNet() [3/3] [DeviceIdNet](#) ([DeviceIdNet](#)% *deviceId*)

10.124.3 Member Function Documentation

10.124.3.1 operator=() [DeviceIdNet](#) operator= ([DeviceIdNet](#)% *deviceId*)

10.124.4 Member Data Documentation

10.124.4.1 BcdDevice `int BcdDevice`

10.124.4.2 BusType `McsBusTypeEnumNet BusType`

10.124.4.3 IdProduct `ProductIdEnumNet IdProduct`

10.124.4.4 IdVendor `VendorIdEnumNet IdVendor`

10.125 DigitalSource< digitalsourceenum > Class Template Reference

Public Member Functions

- [DigitalSource](#) (digitalsourceenum source)
- int [MaxBitNumber](#) ()

Static Public Member Functions

- static int [size](#) ()

Properties

- digitalsourceenum [Source](#) [get]

10.125.1 Constructor & Destructor Documentation

10.125.1.1 DigitalSource() `DigitalSource (digitalsourceenum source)`

10.125.2 Member Function Documentation

10.125.2.1 MaxBitNumber() `int MaxBitNumber ()`

10.125.2.2 size() `static int size () [static]`

10.125.3 Property Documentation

10.125.3.1 Source `digitalsourceenum Source [get]`

10.126 DriverVersionNet Class Reference

Class gives firmware versions of the device's firmware destinations.

Public Member Functions

- [DriverVersionNet](#) ()
Constructor.
- [~DriverVersionNet](#) ()
Destructor.
- unsigned int [GetStatus](#) (CFirmwareDestinationNet dest)
Get status of firmware destination.
- unsigned int [GetStatus](#) (unsigned int index)
Get status of firmware destination.
- unsigned int [GetVersionInt](#) (CFirmwareDestinationNet dest)
Get the version number of firmware destination (major in high word, minor in low word)
- unsigned int [GetVersionInt](#) (unsigned int index)
Get the version number of firmware destination (major in high word, minor in low word)
- unsigned int [GetMajor](#) (CFirmwareDestinationNet dest)
Get the major version number of firmware destination.
- unsigned int [GetMajor](#) (unsigned int index)
Get the major version number of firmware destination.
- unsigned int [GetMinor](#) (CFirmwareDestinationNet dest)
Get the minor version number of firmware destination.
- unsigned int [GetMinor](#) (unsigned int index)
Get the minor version number of firmware destination.
- unsigned int [GetNumEntries](#) ()
Get the number of available firmware destinations.
- String ^ [GetVersionString](#) (CFirmwareDestinationNet dest)
Get the version as a string in the format Major.Minor.
- String ^ [GetVersionString](#) (unsigned int index)
Get the version as a string in the format Major.Minor.
- CFirmwareDestinationNet [GetDestinationCode](#) (unsigned int index)
Get CFirmwareDestinationNet.
- String ^ [GetDestinationName](#) (CFirmwareDestinationNet dest)
Get firmware destination name.
- String ^ [GetDestinationName](#) (unsigned int index)
Get firmware destination name.

Static Public Member Functions

- static String ^ [DriverVersionNet::FormatVersion](#) (unsigned int v)

10.126.1 Detailed Description

Class gives firmware versions of the device's firmware destinations.

10.126.2 Constructor & Destructor Documentation**10.126.2.1 DriverVersionNet()** [DriverVersionNet](#) ()

Constructor.

10.126.2.2 ~DriverVersionNet() [~DriverVersionNet](#) ()

Destructor.

10.126.3 Member Function Documentation**10.126.3.1 DriverVersionNet::FormatVersion()** static String ^ [DriverVersionNet::FormatVersion](#) (unsigned int v) [static]**10.126.3.2 GetDestinationCode()** CFirmwareDestinationNet [GetDestinationCode](#) (unsigned int index)

Get CFirmwareDestinationNet.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.3 GetDestinationName() [1/2] String ^ [GetDestinationName](#) (CFirmwareDestinationNet dest)

Get firmware destination name.

Parameters

<i>dest</i>	by CFirmwareDestiationNet
-------------	---------------------------

10.126.3.4 GetDestinationName() [2/2] `String ^ GetDestinationName (unsigned int index)`

Get firmware destination name.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.5 GetMajor() [1/2] `unsigned int GetMajor (CFirmwareDestiationNet dest)`

Get the major version number of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestiationNet
-------------	---------------------------

10.126.3.6 GetMajor() [2/2] `unsigned int GetMajor (unsigned int index)`

Get the major version number of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.7 GetMinor() [1/2] `unsigned int GetMinor (CFirmwareDestiationNet dest)`

Get the minor version number of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestiationNet
-------------	---------------------------

10.126.3.8 GetMinor() [2/2] unsigned int GetMinor (
 unsigned int *index*)

Get the minor version number of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.9 GetNumEntries() unsigned int GetNumEntries ()

Get the number of available firmware destinations.

10.126.3.10 GetStatus() [1/2] unsigned int GetStatus (
 CFirmwareDestinationNet *dest*)

Get status of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

10.126.3.11 GetStatus() [2/2] unsigned int GetStatus (
 unsigned int *index*)

Get status of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.12 GetVersionInt() [1/2] unsigned int GetVersionInt (
 CFirmwareDestinationNet *dest*)

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

10.126.3.13 GetVersionInt() [2/2] unsigned int GetVersionInt (unsigned int *index*)

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

10.126.3.14 GetVersionString() [1/2] String ^ GetVersionString (CFirmwareDestinationNet *dest*)

Get the version as a string in the format Major.Minor.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

10.126.3.15 GetVersionString() [2/2] String ^ GetVersionString (unsigned int *index*)

Get the version as a string in the format Major.Minor.

Parameters

<i>index</i>	by index of firmware
--------------	----------------------

10.127 FirmwareDestinationNames Class Reference

Static Public Attributes

- static String ^ **DSP** = gcnew String("DSP")
- static String ^ **USB** = gcnew String("USB")
- static String ^ **MCU1** = gcnew String("MCU1")
- static String ^ **Bootstrap** = gcnew String("Bootstrap")
- static String ^ **MCSBUS1** = gcnew String("McsBus1")

- static String ^ MCSBUS2 = gcnew String("McsBus2")
- static String ^ MCSBUS3 = gcnew String("McsBus3")
- static String ^ MCSBUS4 = gcnew String("McsBus4")
- static String ^ MCSBUS5 = gcnew String("McsBus5")
- static String ^ MCSBUS6 = gcnew String("McsBus6")
- static String ^ MCSBUS7 = gcnew String("McsBus7")
- static String ^ MCSBUS8 = gcnew String("McsBus8")
- static String ^ MCSBUS9 = gcnew String("McsBus9")
- static String ^ MCSBUS10 = gcnew String("McsBus10")
- static String ^ MCSBUS11 = gcnew String("McsBus11")
- static String ^ MCSBUS12 = gcnew String("McsBus12")
- static String ^ MCSBUS13 = gcnew String("McsBus13")
- static String ^ BUS1_MCSBUS1 = gcnew String("Bus1McsBus1")
- static String ^ BUS1_MCSBUS2 = gcnew String("Bus1McsBus2")
- static String ^ PIC = gcnew String("PIC")
- static String ^ PIC2 = gcnew String("PIC2")
- static String ^ PIC3 = gcnew String("PIC3")
- static String ^ PIC4 = gcnew String("PIC4")
- static String ^ Altera = gcnew String("Altera")
- static String ^ FPGA2 = gcnew String("FPGA2")
- static String ^ FPGA3 = gcnew String("FPGA3")
- static String ^ FPGA4 = gcnew String("FPGA4")
- static String ^ FPGA5 = gcnew String("FPGA5")
- static String ^ FPGA6 = gcnew String("FPGA6")

10.127.1 Member Data Documentation

10.127.1.1 Altera String ^ Altera = gcnew String("Altera") [static]

10.127.1.2 Bootstrap String ^ Bootstrap = gcnew String("Bootstrap") [static]

10.127.1.3 BUS1_MCSBUS1 String ^ BUS1_MCSBUS1 = gcnew String("Bus1McsBus1") [static]

10.127.1.4 BUS1_MCSBUS2 String ^ BUS1_MCSBUS2 = gcnew String("Bus1McsBus2") [static]

10.127.1.5 DSP String ^ DSP = gcnew String("DSP") [static]

10.127.1.6 FPGA2 `String ^ FPGA2 = gcnew String("FPGA2") [static]`

10.127.1.7 FPGA3 `String ^ FPGA3 = gcnew String("FPGA3") [static]`

10.127.1.8 FPGA4 `String ^ FPGA4 = gcnew String("FPGA4") [static]`

10.127.1.9 FPGA5 `String ^ FPGA5 = gcnew String("FPGA5") [static]`

10.127.1.10 FPGA6 `String ^ FPGA6 = gcnew String("FPGA6") [static]`

10.127.1.11 MCSBUS1 `String ^ MCSBUS1 = gcnew String("McsBus1") [static]`

10.127.1.12 MCSBUS10 `String ^ MCSBUS10 = gcnew String("McsBus10") [static]`

10.127.1.13 MCSBUS11 `String ^ MCSBUS11 = gcnew String("McsBus11") [static]`

10.127.1.14 MCSBUS12 `String ^ MCSBUS12 = gcnew String("McsBus12") [static]`

10.127.1.15 MCSBUS13 `String ^ MCSBUS13 = gcnew String("McsBus13") [static]`

10.127.1.16 MCSBUS2 `String ^ MCSBUS2 = gcnew String("McsBus2") [static]`

10.127.1.17 MCSBUS3 `String ^ MCSBUS3 = gcnew String("McsBus3") [static]`

10.127.1.18 MCSBUS4 `String ^ MCSBUS4 = gcnew String("McsBus4") [static]`

10.127.1.19 MCSBUS5 `String ^ MCSBUS5 = gcnew String("McsBus5") [static]`

10.127.1.20 MCSBUS6 `String ^ MCSBUS6 = gcnew String("McsBus6") [static]`

10.127.1.21 MCSBUS7 `String ^ MCSBUS7 = gcnew String("McsBus7") [static]`

10.127.1.22 MCSBUS8 `String ^ MCSBUS8 = gcnew String("McsBus8") [static]`

10.127.1.23 MCSBUS9 `String ^ MCSBUS9 = gcnew String("McsBus9") [static]`

10.127.1.24 MCU1 `String ^ MCU1 = gcnew String("MCU1") [static]`

10.127.1.25 PIC `String ^ PIC = gcnew String("PIC") [static]`

10.127.1.26 PIC2 `String ^ PIC2 = gcnew String("PIC2") [static]`


10.127.1.27 PIC3 `String ^ PIC3 = gcnew String("PIC3") [static]`

10.127.1.28 PIC4 `String ^ PIC4 = gcnew String("PIC4") [static]`

10.127.1.29 USB `String ^ USB = gcnew String("USB") [static]`

10.128 HeadStageIDType Class Reference

Inheritance diagram for HeadStageIDType:

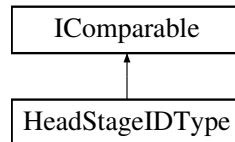

Public Types

- enum [HeadstageTypeEnum](#) {
[Unknown](#),
[MeasuringOnly](#),
[OpticalStimulation](#),
[ElectricalStimulation](#) }

Public Member Functions

- [HeadStageIDType](#) (unsigned int entry, [CW2100_FunctionNet](#)^ device)
- virtual `System::String ^ ToString ()` override
- virtual `bool Equals (Object^ obj)` override
- virtual `Int32 CompareTo (Object^ obj)`

Properties

- `bool Valid` [get]
- `unsigned int Entry` [get]
- `unsigned short ID` [get]
- `System::String^ SN` [get]
- `unsigned int TypeValue` [get]
- `System::String^ Type` [get]
- [HeadstageTypeEnum](#) [HeadstageType](#) [get]
- `System::String^ UserDefinedName` [get]
- `int NumberOfAnalogChannels` [get]
- `int NumberOfStimulationChannels` [get]
- [W2100_StimulusParametersNet](#)^ [StimulusParameters](#) [get]
- `bool HasIMU` [get]
- `bool W16IsW14` [get]
- `bool HasOptoCurrentMeasurement` [get]

10.128.1 Member Enumeration Documentation

10.128.1.1 HeadstageTypeEnum enum [HeadstageTypeEnum](#) [strong]

Enumerator

Unknown	
MeasuringOnly	
OpticalStimulation	
ElectricalStimulation	

10.128.2 Constructor & Destructor Documentation

10.128.2.1 HeadStageIDType() `HeadStageIDType (unsigned int entry, CW2100_FunctionNet^ device)`

10.128.3 Member Function Documentation

10.128.3.1 CompareTo() `virtual Int32 CompareTo (Object^ obj) [virtual]`

10.128.3.2 Equals() `virtual bool Equals (Object^ obj) [override], [virtual]`

10.128.3.3 ToString() `virtual System::String ^ ToString () [override], [virtual]`

10.128.4 Property Documentation

10.128.4.1 Entry `unsigned int Entry [get]`

10.128.4.2 HasIMU `bool HasIMU [get]`

10.128.4.3 HasOptoCurrentMessurement bool HasOptoCurrentMessurement [get]

10.128.4.4 HeadstageType [HeadstageTypeEnum](#) HeadstageType [get]

10.128.4.5 ID unsigned short ID [get]

10.128.4.6 NumberOfAnalogChannels int NumberOfAnalogChannels [get]

10.128.4.7 NumberOfStimulationChannels int NumberOfStimulationChannels [get]

10.128.4.8 SN System:: String^ SN [get]

10.128.4.9 StimulusParameters [W2100_StimulusParametersNet](#)^ StimulusParameters [get]

10.128.4.10 Type System:: String^ Type [get]

10.128.4.11 TypeValue unsigned int TypeValue [get]

10.128.4.12 UserDefinedName System:: String^ UserDefinedName [get]

10.128.4.13 Valid bool Valid [get]

10.128.4.14 W16IsW14 bool W16IsW14 [get]

10.129 HeadstageIDTypeObject Class Reference

Public Member Functions

- [HeadstageIDTypeObject](#) ([HeadStageIDType](#)^ idType)
- virtual [String](#) ^ [ToString](#) () override
- virtual bool [Equals](#) ([Object](#)^ obj) override
- virtual int [GetHashCode](#) () override

Public Attributes

- [HeadStageIDType](#) ^ [_IdType](#)
- [String](#) ^ [_AdditionalText](#)

Properties

- [HeadStageIDType](#)^ [IdType](#) [get]
- [String](#)^ [AdditionalText](#) [get, set]

10.129.1 Constructor & Destructor Documentation

10.129.1.1 HeadstageIDTypeObject() [HeadstageIDTypeObject](#) (
 [HeadStageIDType](#)^ idType)

10.129.2 Member Function Documentation

10.129.2.1 Equals() virtual bool [Equals](#) (
 [Object](#)^ obj) [override], [virtual]

10.129.2.2 GetHashCode() virtual int [GetHashCode](#) () [override], [virtual]

10.129.2.3 ToString() virtual [String](#) ^ [ToString](#) () [override], [virtual]

10.129.3 Member Data Documentation

10.129.3.1 `_AdditionalText` `String ^ _AdditionalText`

10.129.3.2 `_IdType` `HeadStageIDType ^ _IdType`

10.129.4 Property Documentation

10.129.4.1 `AdditionalText` `String^ AdditionalText [get], [set]`

10.129.4.2 `IdType` `HeadStageIDType^ IdType [get]`

10.130 HeadStageIDTypeState Class Reference

Properties

- unsigned int `State` `[get]`
- `HeadStageIDType^ IdType` `[get]`
- bool `ControlState` `[get]`
- bool `DataState` `[get]`

10.130.1 Property Documentation

10.130.1.1 `ControlState` `bool ControlState [get]`

10.130.1.2 `DataState` `bool DataState [get]`

10.130.1.3 `IdType` `HeadStageIDType^ IdType [get]`

10.130.1.4 `State` `unsigned int State [get]`

10.131 mkfilterNet Class Reference

Static Public Member Functions

- static int [mkfilter](#) (String[^] *filtertype*, double *value*, String[^] *passtype*, int *order*, double *alpha1*, double *alpha2*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Amplification*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS_k](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Amplification*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *coeffs*)
- static int [mkfilter_MCS_k](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *coeffs*)
- static void [mkfilter_coef_in_one_set](#) (int *n*, [System::Runtime::InteropServices::In] array< double >[^] *xcoeffs*, [System::Runtime::InteropServices::In] array< double >[^] *ycoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *out_coeffs*)
- static void [mkfilter_scale_coef_in_one_set](#) (int *n*, double *scale*, [System::Runtime::InteropServices::In] array< double >[^] *xcoeffs*, [System::Runtime::InteropServices::In] array< double >[^] *ycoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *out_coeffs*)
- static void [mkfilter_normalize_coeffs_short](#) (short *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< short >[^]% *out_coeffs*)
- static void [mkfilter_normalize_coeffs_int](#) (int *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< int >[^]% *out_coeffs*)
- static void [mkfilter_normalize_scale_coeffs_int](#) (int *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< int >[^]% *out_coeffs*)
- static double [mkfilter_highpass_coeff](#) (int *SamplesPerSecond*, double *Frequency*)
- static double [mkfilter_highpass_k](#) (int *SamplesPerSecond*, double *Frequency*)
- static double [mkfilter_highpass_frequency_from_coeff](#) (int *SamplesPerSecond*, double *coeff*)
- static double [mkfilter_highpass_frequency_from_k](#) (int *SamplesPerSecond*, double *k*)

10.131.1 Member Function Documentation

10.131.1.1 [mkfilter\(\)](#) static int [mkfilter](#) (
String[^] *filtertype*,
double *value*,
String[^] *passtype*,
int *order*,
double *alpha1*,
double *alpha2*,
[System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*,
[System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*) [static]

10.131.1.2 mkfilter_coef_in_one_set() static void mkfilter_coef_in_one_set (
 int *n*,
 [System::Runtime::InteropServices::In] array< double >^ *xcoeffs*,
 [System::Runtime::InteropServices::In] array< double >^ *ycoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *out_coeffs*) [static]

10.131.1.3 mkfilter_highpass_coeff() static double mkfilter_highpass_coeff (
 int *SamplesPerSecond*,
 double *Frequency*) [static]

10.131.1.4 mkfilter_highpass_frequency_from_coeff() static double mkfilter_highpass_frequency_from_↵
 from_coeff (
 int *SamplesPerSecond*,
 double *coeff*) [static]

10.131.1.5 mkfilter_highpass_frequency_from_k() static double mkfilter_highpass_frequency_from_↵
 _k (
 int *SamplesPerSecond*,
 double *k*) [static]

10.131.1.6 mkfilter_highpass_k() static double mkfilter_highpass_k (
 int *SamplesPerSecond*,
 double *Frequency*) [static]

10.131.1.7 mkfilter_MCS() [1/2] static int mkfilter_MCS (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Amplification*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *xcoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *ycoeffs*) [static]

10.131.1.8 mkfilter_MCS() [2/2] static int mkfilter_MCS (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *xcoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *ycoeffs*) [static]

- 10.131.1.9 mkfilter_MCS_k()** [1/2] static int mkfilter_MCS_k (
int *SamplesPerSecond*,
double *R1*,
double *R2*,
double *C*,
double *Amplification*,
double *Correction*,
[System::Runtime::InteropServices::Out] array< double >^% *coeffs*) [static]
- 10.131.1.10 mkfilter_MCS_k()** [2/2] static int mkfilter_MCS_k (
int *SamplesPerSecond*,
double *R1*,
double *R2*,
double *C*,
double *Correction*,
[System::Runtime::InteropServices::Out] array< double >^% *coeffs*) [static]
- 10.131.1.11 mkfilter_normalize_coeffs_int()** static void mkfilter_normalize_coeffs_int (
int *maxvalue*,
[System::Runtime::InteropServices::In] array< double >^ *coeffs*,
[System::Runtime::InteropServices::Out] array< int >^% *out_coeffs*) [static]
- 10.131.1.12 mkfilter_normalize_coeffs_short()** static void mkfilter_normalize_coeffs_short (
short *maxvalue*,
[System::Runtime::InteropServices::In] array< double >^ *coeffs*,
[System::Runtime::InteropServices::Out] array< short >^% *out_coeffs*) [static]
- 10.131.1.13 mkfilter_normalize_scale_coeffs_int()** static void mkfilter_normalize_scale_coeffs_int (
int *maxvalue*,
[System::Runtime::InteropServices::In] array< double >^ *coeffs*,
[System::Runtime::InteropServices::Out] array< int >^% *out_coeffs*) [static]
- 10.131.1.14 mkfilter_scale_coef_in_one_set()** static void mkfilter_scale_coef_in_one_set (
int *n*,
double *scale*,
[System::Runtime::InteropServices::In] array< double >^ *xcoeffs*,
[System::Runtime::InteropServices::In] array< double >^ *ycoeffs*,
[System::Runtime::InteropServices::Out] array< double >^% *out_coeffs*) [static]

10.132 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference

Public Member Functions

- void [SetParameter](#) (unsigned short command, unsigned short index, unsigned int value)
- void [SetParameter](#) (unsigned short command, unsigned short index, unsigned int value1, unsigned int value2)
- void [SetUserParameter](#) (unsigned short index, unsigned int value)
 - Stores persistently 32 bit integer values on motor controller*
- void [SetUserParameter](#) (unsigned short index, int value)
 - Stores persistently 32 bit integer values on motor controller*
- void [GetParameter](#) (unsigned short command, unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value)
- void [GetParameter](#) (unsigned short command, unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value1, [System::Runtime::InteropServices::Out]unsigned int% value2)
- void [GetUserParameter](#) (unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value)
 - Reads 32 bit integer values stored persistently on motor controller*
- void [GetUserParameter](#) (unsigned short index, [System::Runtime::InteropServices::Out]int% value)
 - Reads 32 bit integer values stored persistently on motor controller*
- void [FindReferencePhase0](#) (unsigned char busaddress, char axes)
- void [FindReferencePhase0](#) (unsigned char busaddress, char axes, int timeout)
- unsigned char [HasRef](#) (unsigned char busaddress, char axes)
- void [SetHWRevision](#) (unsigned int revision)
- unsigned int [GetHWRevision](#) ()
- void [SetHWConfig](#) (unsigned int config)
- unsigned int [GetHWConfig](#) ()
- void [SetMinPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMinPressureWaitTime](#) ()
- void [SetMinPressure](#) (unsigned int pressure)
- unsigned int [GetMinPressure](#) ()
- void [SetMaxPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMaxPressureWaitTime](#) ()
- void [SetMinNoPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMinNoPressureWaitTime](#) ()
- void [SetMaxNoPressure](#) (unsigned int pressure)
- unsigned int [GetMaxNoPressure](#) ()
- void [SetMaxNoPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMaxNoPressureWaitTime](#) ()
- void [SetSearchReferenceMethod](#) (unsigned char busaddress, char axes, unsigned int method)
- unsigned int [GetSearchReferenceMethod](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceOffsetPos](#) (unsigned char busaddress, char axes, int offsetpos)
- int [GetSearchReferenceOffsetPos](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFastSpeed](#) (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short [GetSearchReferenceFastSpeed](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFastAccel](#) (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short [GetSearchReferenceFastAccel](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFineSpeed](#) (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short [GetSearchReferenceFineSpeed](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFineAccel](#) (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short [GetSearchReferenceFineAccel](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceMoveOut](#) (unsigned char busaddress, char axes, int move)
- int [GetSearchReferenceMoveOut](#) (unsigned char busaddress, char axes)
- void [SetAxisConfig](#) (unsigned char busaddress, char axes, unsigned int config)
- unsigned int [GetAxisConfig](#) (unsigned char busaddress, char axes)
- void [GetPhases](#) (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out] unsigned short% phase0, [System::Runtime::InteropServices::Out] unsigned short% lastphase)

10.132.1 Member Function Documentation

10.132.1.1 FindReferencePhase0() [1/2] void FindReferencePhase0 (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.2 FindReferencePhase0() [2/2] void FindReferencePhase0 (
 unsigned char *busaddress*,
 char *axes*,
 int *timeout*)

10.132.1.3 GetAxisConfig() unsigned int GetAxisConfig (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.4 GetHWConfig() unsigned int GetHWConfig ()

10.132.1.5 GetHWRevision() unsigned int GetHWRevision ()

10.132.1.6 GetMaxNoPressure() unsigned int GetMaxNoPressure ()

10.132.1.7 GetMaxNoPressureWaitTime() unsigned int GetMaxNoPressureWaitTime ()

10.132.1.8 GetMaxPressureWaitTime() unsigned int GetMaxPressureWaitTime ()

10.132.1.9 GetMinNoPressureWaitTime() unsigned int GetMinNoPressureWaitTime ()

10.132.1.10 GetMinPressure() unsigned int GetMinPressure ()

10.132.1.11 GetMinPressureWaitTime() unsigned int GetMinPressureWaitTime ()

10.132.1.12 GetParameter() [1/2] void GetParameter (
 unsigned short *command*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned int% *value*)

10.132.1.13 GetParameter() [2/2] void GetParameter (
 unsigned short *command*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned int% *value1*,
 [System::Runtime::InteropServices::Out] unsigned int% *value2*)

10.132.1.14 GetPhases() void GetPhases (
 unsigned char *busaddress*,
 char *axes*,
 [System::Runtime::InteropServices::Out] unsigned short% *phase0*,
 [System::Runtime::InteropServices::Out] unsigned short% *lastphase*)

10.132.1.15 GetSearchReferenceFastAccel() unsigned short GetSearchReferenceFastAccel (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.16 GetSearchReferenceFastSpeed() unsigned short GetSearchReferenceFastSpeed (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.17 GetSearchReferenceFineAccel() unsigned short GetSearchReferenceFineAccel (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.18 GetSearchReferenceFineSpeed() unsigned short GetSearchReferenceFineSpeed (unsigned char *busaddress*, char *axes*)

10.132.1.19 GetSearchReferenceMethod() unsigned int GetSearchReferenceMethod (unsigned char *busaddress*, char *axes*)

10.132.1.20 GetSearchReferenceMoveOut() int GetSearchReferenceMoveOut (unsigned char *busaddress*, char *axes*)

10.132.1.21 GetSearchReferenceOffsetPos() int GetSearchReferenceOffsetPos (unsigned char *busaddress*, char *axes*)

10.132.1.22 GetUserParameter() [1/2] void GetUserParameter (unsigned short *index*, [System::Runtime::InteropServices::Out] int% *value*)

Reads 32 bit integer values stored persistently on motor controller

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data buffer

10.132.1.23 GetUserParameter() [2/2] void GetUserParameter (unsigned short *index*, [System::Runtime::InteropServices::Out] unsigned int% *value*)

Reads 32 bit integer values stored persistently on motor controller

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data buffer

10.132.1.24 HasRef() unsigned char HasRef (
 unsigned char *busaddress*,
 char *axes*)

10.132.1.25 SetAxisConfig() void SetAxisConfig (
 unsigned char *busaddress*,
 char *axes*,
 unsigned int *config*)

10.132.1.26 SetHWConfig() void SetHWConfig (
 unsigned int *config*)

10.132.1.27 SetHWRevision() void SetHWRevision (
 unsigned int *revision*)

10.132.1.28 SetMaxNoPressure() void SetMaxNoPressure (
 unsigned int *pressure*)

10.132.1.29 SetMaxNoPressureWaitTime() void SetMaxNoPressureWaitTime (
 unsigned int *t*)

10.132.1.30 SetMaxPressureWaitTime() void SetMaxPressureWaitTime (
 unsigned int *t*)

10.132.1.31 SetMinNoPressureWaitTime() void SetMinNoPressureWaitTime (
 unsigned int *t*)

10.132.1.32 SetMinPressure() void SetMinPressure (
 unsigned int *pressure*)

10.132.1.33 SetMinPressureWaitTime() `void SetMinPressureWaitTime (unsigned int t)`

10.132.1.34 SetParameter() [1/2] `void SetParameter (unsigned short command, unsigned short index, unsigned int value)`

10.132.1.35 SetParameter() [2/2] `void SetParameter (unsigned short command, unsigned short index, unsigned int value1, unsigned int value2)`

10.132.1.36 SetSearchReferenceFastAccel() `void SetSearchReferenceFastAccel (unsigned char busaddress, char axes, unsigned short accel)`

10.132.1.37 SetSearchReferenceFastSpeed() `void SetSearchReferenceFastSpeed (unsigned char busaddress, char axes, unsigned short speed)`

10.132.1.38 SetSearchReferenceFineAccel() `void SetSearchReferenceFineAccel (unsigned char busaddress, char axes, unsigned short accel)`

10.132.1.39 SetSearchReferenceFineSpeed() `void SetSearchReferenceFineSpeed (unsigned char busaddress, char axes, unsigned short speed)`

10.132.1.40 SetSearchReferenceMethod() `void SetSearchReferenceMethod (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned int method)`

10.132.1.41 SetSearchReferenceMoveOut() `void SetSearchReferenceMoveOut (`
 `unsigned char busaddress,`
 `char axes,`
 `int move)`

10.132.1.42 SetSearchReferenceOffsetPos() `void SetSearchReferenceOffsetPos (`
 `unsigned char busaddress,`
 `char axes,`
 `int offsetpos)`

10.132.1.43 SetUserParameter() [1/2] `void SetUserParameter (`
 `unsigned short index,`
 `int value)`

Stores *persistently* 32 bit integer values on motor controller

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data to be stored

10.132.1.44 SetUserParameter() [2/2] `void SetUserParameter (`
 `unsigned short index,`
 `unsigned int value)`

Stores *persistently* 32 bit integer values on motor controller

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data to be stored

10.133 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference

Public Member Functions

- void [FindReferencePhase0XY](#) ()
- void [FindReferencePhase0XY](#) (int timeout)

10.133.1 Member Function Documentation

10.133.1.1 [FindReferencePhase0XY\(\)](#) [1/2] `void FindReferencePhase0XY ()`

10.133.1.2 [FindReferencePhase0XY\(\)](#) [2/2] `void FindReferencePhase0XY (
int timeout)`

10.134 CMeaAudioFunctionNet::s_setaudionet Struct Reference

Public Attributes

- int [channel](#)
- int [amplification](#)

10.134.1 Member Data Documentation

10.134.1.1 [amplification](#) `int amplification`

10.134.1.2 [channel](#) `int channel`

10.135 CStimulusFunctionNet::SidebandData Class Reference

Public Member Functions

- [SidebandData](#) ()
- [~SidebandData](#) ()
Destructor: called by Dispose()
- [!SidebandData](#) ()
Finalizer: called by GC before collecting

Properties

- `array< int32_t >^ Sideband` [get]
- `array< uint64_t >^ Duration` [get]

10.135.1 Constructor & Destructor Documentation

10.135.1.1 SidebandData() `SidebandData ()`

10.135.1.2 ~SidebandData() `~SidebandData ()`

Destructor: called by Dispose()

10.135.1.3 "!SidebandData() `!SidebandData ()`

Finalizer: called by GC before collecting

10.135.2 Property Documentation

10.135.2.1 Duration `array< uint64_t >^ Duration` [get]

10.135.2.2 Sideband `array< int32_t >^ Sideband` [get]

10.136 StgStatusNet Class Reference

Static Public Member Functions

- static `StgStatusNet ^ FromIntPtr` (IntPtr stgstatus)
- static `StgStatusNet ^ FromPtr` (stgstatus_t *stgstatus)

Public Attributes

- `array< Stg200xTriggerStatusEnumNet > ^ TiggerStatus`
- `array< uint32_t > ^ ListOfChangedTriggers`

10.136.1 Member Function Documentation

10.136.1.1 FromIntPtr() static [StgStatusNet](#) ^ FromIntPtr (
 IntPtr stgstatus) [static]

10.136.1.2 FromPtr() static [StgStatusNet](#) ^ FromPtr (
 stgstatus_t * stgstatus) [static]

10.136.2 Member Data Documentation

10.136.2.1 ListOfChangedTriggers array<uint32_t> ^ ListOfChangedTriggers

10.136.2.2 TiggerStatus array<Stg200xTriggerStatusEnumNet> ^ TiggerStatus

10.137 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference

Public Member Functions

- [StimulusDeviceDataAndUnrolledData](#) ()
- [~StimulusDeviceDataAndUnrolledData](#) ()
Destructor: called by Dispose()
- [!StimulusDeviceDataAndUnrolledData](#) ()
Finalizer: called by GC before collecting

Properties

- array< uint8_t >^ [DeviceData](#) [get]
- int [DeviceDataLength](#) [get]
- array< int32_t >^ [UnrolledAmplitude](#) [get]
- array< uint32_t >^ [UnrolledSync](#) [get]
- array< uint64_t >^ [UnrolledDuration](#) [get]

10.137.1 Constructor & Destructor Documentation

10.137.1.1 StimulusDeviceDataAndUnrolledData() [StimulusDeviceDataAndUnrolledData \(\)](#)

10.137.1.2 ~StimulusDeviceDataAndUnrolledData() [~StimulusDeviceDataAndUnrolledData \(\)](#)

Destructor: called by Dispose()

10.137.1.3 ~!StimulusDeviceDataAndUnrolledData() [!StimulusDeviceDataAndUnrolledData \(\)](#)

Finalizer: called by GC before collecting

10.137.2 Property Documentation

10.137.2.1 DeviceData [array< uint8_t>^ DeviceData \[get\]](#)

10.137.2.2 DeviceDataLength [int DeviceDataLength \[get\]](#)

10.137.2.3 UnrolledAmplitude [array< int32_t>^ UnrolledAmplitude \[get\]](#)

10.137.2.4 UnrolledDuration [array< uint64_t>^ UnrolledDuration \[get\]](#)

10.137.2.5 UnrolledSync [array< uint32_t>^ UnrolledSync \[get\]](#)

10.138 usbSetupPacket_t Class Reference

Public Attributes

- [uint8_t bmRequestType](#)
- [uint8_t bRequest](#)
- [uint16_t wValue](#)
- [uint16_t wIndex](#)
- [uint16_t wLength](#)

10.138.1 Member Data Documentation

10.138.1.1 bmRequestType `uint8_t bmRequestType`

10.138.1.2 bRequest `uint8_t bRequest`

10.138.1.3 wIndex `uint16_t wIndex`

10.138.1.4 wLength `uint16_t wLength`

10.138.1.5 wValue `uint16_t wValue`

10.139 W2100_StimulusParametersNet Struct Reference

Public Attributes

- `int DACResolution`
- `int TimeResolutionInNanoSeconds`
- `int VoltageRangeInMicroVolt`
- `int VoltageResolutionInMicroVolt`
- `int CurrentRangeInNanoAmp`
- `int CurrentResolutionInNanoAmp`

10.139.1 Member Data Documentation

10.139.1.1 CurrentRangeInNanoAmp `int CurrentRangeInNanoAmp`

10.139.1.2 CurrentResolutionInNanoAmp `int CurrentResolutionInNanoAmp`

10.139.1.3 DACResolution `int DACResolution`

10.139.1.4 TimeResolutionInNanoSeconds `int TimeResolutionInNanoSeconds`

10.139.1.5 VoltageRangeInMicroVolt `int VoltageRangeInMicroVolt`

10.139.1.6 VoltageResolutionInMicroVolt `int VoltageResolutionInMicroVolt`

Index

- !CDacCalibrationFunctionNet
 - CDacCalibrationFunctionNet, [48](#)
- !CDigOutStimulatorFunctionNet
 - CDigOutStimulatorFunctionNet, [58](#)
- !CExternDTesterDeviceNet
 - CExternDTesterDeviceNet, [62](#)
- !CInterfaceboardFunctionNet
 - CInterfaceboardFunctionNet, [106](#)
- !CLIH3DeviceNet
 - CLIH3DeviceNet, [109](#)
- !CMEA2100x256FunctionNet
 - CMEA2100x256FunctionNet, [246](#)
- !CMcsUsbFunctionNet
 - CMcsUsbFunctionNet, [216](#)
- !CMcsUsbListNet
 - CMcsUsbListNet, [222](#)
- !CMcsUsbNet
 - CMcsUsbNet, [228](#)
- !CMeFunctionNet
 - CMeFunctionNet, [280](#)
- !CMeaCleanDeviceNet
 - CMeaCleanDeviceNet, [252](#)
- !CMeaCoatDeviceNet
 - CMeaCoatDeviceNet, [256](#)
- !CMultiBatteryChargerDeviceNet
 - CMultiBatteryChargerDeviceNet, [282](#)
- !CMultiwellCallbackFunctionNet
 - CMultiwellCallbackFunctionNet, [289](#)
- !CMultiwellDeviceNet
 - CMultiwellDeviceNet, [293](#)
- !CMultiwellOptoStimFunctionNet
 - CMultiwellOptoStimFunctionNet, [298](#)
- !CPPPCFunctionNet
 - CPPPCFunctionNet, [325](#)
- !CPedoterDeviceNet
 - CPedoterDeviceNet, [312](#)
- !CPositionIIDeviceNet
 - CPositionIIDeviceNet, [317](#)
- !CPositionImpDeviceNet
 - CPositionImpDeviceNet, [320](#)
- !CProgramPressureCurveNet
 - CProgramPressureCurveNet, [337](#)
- !CPulseGeneratorFunctionNet
 - CPulseGeneratorFunctionNet, [339](#)
- !CRFFFunctionNet
 - CRFFFunctionNet, [347](#)
- !CSCUFunctionNet
 - CSCUFunctionNet, [399](#)
- !CTEERFunctionNet
 - CTEERFunctionNet, [489](#)
- !CWarnerUssingDeviceNet
 - CWarnerUssingDeviceNet, [515](#)
- !CWarnerUssingFunctionNet
 - CWarnerUssingFunctionNet, [518](#)
- !CWarnerValveControllerDeviceNet
 - CWarnerValveControllerDeviceNet, [536](#)
- !CWarnerValveControllerDeviceTesterFunctionNet
 - CWarnerValveControllerDeviceTesterFunctionNet, [556](#)
- !SidebandData
 - CStimulusFunctionNet::SidebandData, [591](#)
- !StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- __AdditionalText
 - HeadstageIDTypeObject, [578](#)
- __IdType
 - HeadstageIDTypeObject, [579](#)
- ~CCMOSMeaDeviceNet
 - CCMOSMeaDeviceNet, [44](#)
- ~CChannelTestDeviceNet
 - CChannelTestDeviceNet, [31](#)
- ~CCreateFilterNet
 - CCreateFilterNet, [46](#)
- ~CDacCalibrationFunctionNet
 - CDacCalibrationFunctionNet, [48](#)
- ~CDigOutStimulatorFunctionNet
 - CDigOutStimulatorFunctionNet, [57](#)
- ~CExternDTesterDeviceNet
 - CExternDTesterDeviceNet, [62](#)
- ~CFilterCoefficientsNet
 - CFilterCoefficientsNet, [65](#)
- ~CFilterPropertyNet
 - CFilterPropertyNet, [70](#)
- ~CFluidControlDeviceNet
 - CFluidControlDeviceNet, [73](#)
- ~CGenericDevelopDeviceNet
 - CGenericDevelopDeviceNet, [86](#)
- ~CGilsonDeviceNet
 - CGilsonDeviceNet, [97](#)
- ~CInterfaceboardFunctionNet
 - CInterfaceboardFunctionNet, [106](#)
- ~CLIH3DeviceNet
 - CLIH3DeviceNet, [109](#)
- ~CMEA2100x256FunctionNet
 - CMEA2100x256FunctionNet, [246](#)
- ~CMcsBusNet
 - CMcsBusNet, [153](#)
- ~CMcsBus_AxisParametersNet
 - CMcsBus_AxisParametersNet, [115](#)
- ~CMcsBus_ExtensionNet
 - CMcsBus_ExtensionNet, [116](#)
- ~CMcsBus_FYIExtensionNet
 - CMcsBus_FYIExtensionNet, [117](#)
- ~CMcsBus_MotorControlNet
 - CMcsBus_MotorControlNet, [122](#)
- ~CMcsBus_SensorNet
 - CMcsBus_SensorNet, [138](#)
- ~CMcsBus_TempSensorNet
 - CMcsBus_TempSensorNet, [147](#)

- ~CMcsBus_VoltageModeNet
 - CMcsBus_VoltageModeNet, [149](#)
- ~CMcsUsbDacqNet
 - CMcsUsbDacqNet, [163](#)
- ~CMcsUsbFactoryNet
 - CMcsUsbFactoryNet, [209](#)
- ~CMcsUsbFunctionNet
 - CMcsUsbFunctionNet, [216](#)
- ~CMcsUsbListEntryNet
 - CMcsUsbListEntryNet, [218](#)
- ~CMcsUsbListNet
 - CMcsUsbListNet, [222](#)
- ~CMcsUsbNet
 - CMcsUsbNet, [228](#)
- ~CMeFunctionNet
 - CMeFunctionNet, [280](#)
- ~CMeaCleanDeviceNet
 - CMeaCleanDeviceNet, [251](#)
- ~CMeaCoatDeviceNet
 - CMeaCoatDeviceNet, [256](#)
- ~CMeaDeviceNet
 - CMeaDeviceNet, [262](#)
- ~CMealImpedanceDeviceNet
 - CMealImpedanceDeviceNet, [274](#)
- ~CMeaSwitchDeviceNet
 - CMeaSwitchDeviceNet, [277](#)
- ~CMeaUSBDeviceNet
 - CMeaUSBDeviceNet, [279](#)
- ~CMultiBatteryChargerDeviceNet
 - CMultiBatteryChargerDeviceNet, [282](#)
- ~CMultiwellCallbackFunctionNet
 - CMultiwellCallbackFunctionNet, [289](#)
- ~CMultiwellDeviceNet
 - CMultiwellDeviceNet, [293](#)
- ~CMultiwellOptoStimFunctionNet
 - CMultiwellOptoStimFunctionNet, [298](#)
- ~CNF_GenDeviceNet
 - CNF_GenDeviceNet, [302](#)
- ~CokuvisionStimulatorDeviceNet
 - CokuvisionStimulatorDeviceNet, [307](#)
- ~CPPCFunctionNet
 - CPPCFunctionNet, [325](#)
- ~CPathIdentDeviceNet
 - CPathIdentDeviceNet, [311](#)
- ~CPedoterDeviceNet
 - CPedoterDeviceNet, [312](#)
- ~CPeristalticPumpDeviceNet
 - CPeristalticPumpDeviceNet, [314](#)
- ~CPgaDeviceNet
 - CPgaDeviceNet, [315](#)
- ~CPositionIIDeviceNet
 - CPositionIIDeviceNet, [317](#)
- ~CPositionImpDeviceNet
 - CPositionImpDeviceNet, [320](#)
- ~CProgramPressureCurveNet
 - CProgramPressureCurveNet, [337](#)
- ~CPulseGeneratorFunctionNet
 - CPulseGeneratorFunctionNet, [339](#)
- ~CRFFFunctionNet
 - CRFFFunctionNet, [347](#)
- ~CRetinaLedDeviceNet
 - CRetinaLedDeviceNet, [344](#)
- ~CRoboDeviceNet
 - CRoboDeviceNet, [370](#)
- ~CRoboFluidDeviceNet
 - CRoboFluidDeviceNet, [382](#)
- ~CSCUFunctionNet
 - CSCUFunctionNet, [399](#)
- ~CSafeISDeviceNet
 - CSafeISDeviceNet, [394](#)
- ~CStg200xBasicNet
 - CStg200xBasicNet, [418](#)
- ~CStg200xDownloadNet
 - CStg200xDownloadNet, [454](#)
- ~CSw2to64DeviceNet
 - CSw2to64DeviceNet, [472](#)
- ~CTEERFunctionNet
 - CTEERFunctionNet, [489](#)
- ~CTEERMachineDeviceNet
 - CTEERMachineDeviceNet, [497](#)
- ~CTcxDeviceNet
 - CTcxDeviceNet, [476](#)
- ~CWarnerUssingDeviceNet
 - CWarnerUssingDeviceNet, [515](#)
- ~CWarnerUssingFunctionNet
 - CWarnerUssingFunctionNet, [518](#)
- ~CWarnerValveControllerDeviceNet
 - CWarnerValveControllerDeviceNet, [536](#)
- ~CWarnerValveControllerDeviceTesterFunctionNet
 - CWarnerValveControllerDeviceTesterFunctionNet, [556](#)
- ~DriverVersionNet
 - DriverVersionNet, [567](#)
- ~SidebandData
 - CStimulusFunctionNet::SidebandData, [591](#)
- ~StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- A
 - CFilterCoefficientsNet, [65](#)
- A1
 - CFilterCoefficientsNet, [65](#)
- A2
 - CFilterCoefficientsNet, [66](#)
- AdditionalText
 - HeadstageIDTypeObject, [579](#)
- AddLoopEntry
 - CRetinaLedDeviceNet, [344](#)
- AddSelectedChannelsQueue
 - CMcsUsbDacqNet, [163–165](#)
- AddSoftwareKey
 - CMcsUsbNet, [228](#)
- AddTableEntry
 - CRetinaLedDeviceNet, [344](#)
- Altera
 - FirmwareDestinationNames, [572](#)

- amplification
 - CMeaAudioFunctionNet::s_setaudionet, [590](#)
 - CW2100_FunctionNet::AudioChannelsNet, [27](#)
- AmplifierSettle
 - CIntanMea_FunctionNet, [104](#)
- AnalogGain
 - CMeaDeviceNet, [267](#)
- ApplyGains
 - CPgaDeviceNet, [315](#)
- AreTransistorVoltagesSet
 - CCMOSMea_FunctionNet, [34](#)
- AssociateToThis
 - CMcsUsbNet, [228](#)
- AutomaticAnalogOut
 - CSCUFunctionNet, [399](#)
- Axes_I
 - CRoboDeviceNet, [377](#)
- Axes_X
 - CRoboDeviceNet, [377](#)
- Axes_Y
 - CRoboDeviceNet, [377](#)
- Axes_Z
 - CRoboDeviceNet, [377](#)
- Axis_I
 - CRoboDeviceNet, [377](#)
- Axis_X
 - CRoboDeviceNet, [377](#)
- Axis_Y
 - CRoboDeviceNet, [377](#)
- Axis_Z
 - CRoboDeviceNet, [378](#)
- B
 - CFilterCoefficientsNet, [66](#)
- B0
 - CFilterCoefficientsNet, [66](#)
- B1
 - CFilterCoefficientsNet, [66](#)
- B2
 - CFilterCoefficientsNet, [66](#)
- BatteryState, [28](#)
 - Charge, [28](#)
 - ChargeRegionString, [28](#)
 - ChargeString, [28](#)
 - Voltage, [28](#)
 - VoltageString, [28](#)
- BcdDevice
 - DeviceldNet, [565](#)
- BeginImpedanceCheck
 - CIntanMea_FunctionNet, [104](#)
- BesselFilterHighPassNet, [28](#)
 - BesselFilterHighPassNet, [29](#)
- BesselFilterLowPassNet, [29](#)
 - BesselFilterLowPassNet, [29](#)
- bmRequestType
 - usbSetupPacket_t, [594](#)
- BOOST_BIT
 - CW2100_StimulatorFunctionNet, [513](#)
- Bootstrap
 - FirmwareDestinationNames, [572](#)
- bRequest
 - usbSetupPacket_t, [594](#)
- BurnAdcOffset
 - COctoPotDeviceNet, [303](#)
- BurnDacOffset
 - CDacCalibrationFunctionNet, [49](#)
 - COctoPotDeviceNet, [304](#)
- BUS1_MCSBUS1
 - FirmwareDestinationNames, [572](#)
- BUS1_MCSBUS2
 - FirmwareDestinationNames, [572](#)
- BusType
 - DeviceldNet, [565](#)
- ButterworthFilterHighPassNet, [29](#)
 - ButterworthFilterHighPassNet, [30](#)
- ButterworthFilterLowPassNet, [30](#)
 - ButterworthFilterLowPassNet, [30](#)
- CalibrateThermocouple
 - CFluidControlDeviceNet, [73](#)
 - CTcxDeviceNet, [476](#)
- CancelInternalCalibration
 - CTEERFunctionNet, [489](#)
- CancelPoolLoop
 - CRoboDeviceNet, [371](#)
- CancelPoolLoopAndStopMovement
 - CRoboDeviceNet, [371](#)
- CancelTableLoop
 - CRoboDacqNet, [356](#)
- CancelTableLoopAndStopTable
 - CRoboDacqNet, [356](#)
- CapacityTest
 - CMultiBatteryChargerDeviceNet, [282](#)
- CatchAmpGetAdcMean
 - CMcsBus_SensorNet, [138](#)
- CatchAmpGetAdcValue
 - CMcsBus_SensorNet, [138](#)
- CatchAmpGetAdcValueH
 - CMcsBus_SensorNet, [138](#)
- CatchAmpGetAdcValueL
 - CMcsBus_SensorNet, [138](#)
- CatchAmpGetDacAmplitude
 - CMcsBus_SensorNet, [139](#)
- CatchAmpGetDacEnable
 - CMcsBus_SensorNet, [139](#)
- CatchAmpGetDacOffset
 - CMcsBus_SensorNet, [139](#)
- CatchAmpGetPwmEnable
 - CMcsBus_SensorNet, [139](#)
- CatchAmpSetDacAmplitude
 - CMcsBus_SensorNet, [139](#)
- CatchAmpSetDacEnable
 - CMcsBus_SensorNet, [139](#)
- CatchAmpSetDacOffset
 - CMcsBus_SensorNet, [139](#)
- CatchAmpSetPwmEnable
 - CMcsBus_SensorNet, [139](#)
- CChannelTestDeviceNet, [31](#)

- ~CChannelTestDeviceNet, [31](#)
- CChannelTestDeviceNet, [31](#)
- SetAmplitude, [31](#)
- SetAttenuation, [31](#)
- SetFrequency, [31](#)
- SetWaveform, [32](#)
- CCMOSMea_FunctionNet, [32](#)
 - AreTransistorVoltagesSet, [34](#)
 - CCMOSMea_FunctionNet, [34](#)
 - ClearSTGOutput, [34](#)
 - DetectChipType, [34](#)
 - EnableChannelsInGroup, [34](#), [35](#)
 - GetADCInputOffset, [35](#)
 - GetBath, [35](#)
 - GetBathMode, [35](#)
 - GetEnabledChannelsInGroup, [35](#)
 - GetGate, [35](#)
 - GetGNDI, [35](#)
 - GetGroupADCBits, [35](#)
 - GetGroupChannelBitmaskBySelect, [36](#)
 - GetGroupChannelBitmaskHS1NCBathCurrent, [36](#)
 - GetGroupChannelBitmaskHS1NCCol2Current, [36](#)
 - GetGroupChannelBitmaskHS1NChipTemp, [36](#)
 - GetGroupChannelBitmaskHS1Sidebands, [37](#)
 - GetGroupChannelBitmaskHS1TriggerStatus, [37](#)
 - GetGroupChannelBitmaskIFDigChannels, [37](#)
 - GetGroupChannelBitmaskInterfaceADC, [37](#)
 - GetGroupChannelBitmaskPacketFrameContext, [38](#)
 - GetGroupChannelBitmaskSTG1DACSignal, [38](#)
 - GetGroupDCOffset, [38](#)
 - GetGroupID, [38](#)
 - GetGroupNumberOfChannels, [39](#)
 - GetGroupResolutionPerDigit, [39](#)
 - GetGroupSampleSize, [39](#)
 - GetGroupType, [39](#)
 - GetGroupUnit, [39](#), [40](#)
 - GetMaxNumOfColumns, [40](#)
 - GetNeurochipMemoryData, [40](#)
 - GetNeurochipMemorySize, [40](#)
 - GetNumberOfSupportedGroups, [40](#)
 - GetSourceBulk, [40](#)
 - GetSourceDrain, [40](#)
 - GetSourceGate, [40](#)
 - GetStimulusSites, [41](#)
 - GetVDD3I, [41](#)
 - GetVDDI, [41](#)
 - IsChipPowered, [41](#)
 - IsGateFloating, [41](#)
 - PowerChip, [41](#)
 - SetADCInputOffset, [41](#)
 - SetBath, [41](#)
 - SetBathMode, [41](#)
 - SetGate, [41](#)
 - SetGateFloating, [42](#)
 - SetGateToVOP, [42](#)
 - SetNeurochipMemoryData, [42](#)
 - SetSourceBulk, [42](#)
 - SetSourceDrain, [42](#)
 - SetSourceGate, [42](#)
 - SetStimulusSites, [42](#)
 - UpdateTransistorVoltages, [42](#)
 - VOPSTimerSetResetTimes, [42](#), [43](#)
- CCMOSMeaDeviceNet, [43](#)
 - ~CCMOSMeaDeviceNet, [44](#)
 - CCMOSMeaDeviceNet, [44](#)
 - CMosMea, [45](#)
 - GetAvailableBaseSamplerates, [44](#)
 - GetBaseSamplerate, [44](#)
 - GetChannelDataI16, [44](#)
 - GetChannelDataI32, [44](#)
 - GetChannelDataUI16, [44](#)
 - GetChannelDataUI32, [45](#)
 - GetCMOSDataDictionary, [45](#)
 - SetBaseSamplerate, [45](#)
 - SetRegionOfInterests, [45](#)
 - Stimulus, [45](#)
 - UpdateChannelBlock, [45](#)
- CCMOSMeaDeviceNet::CRegionOfInterestRect, [343](#)
 - CRegionOfInterestRect, [343](#)
 - DeepCopy, [343](#)
 - m_Bottom, [343](#)
 - m_Left, [343](#)
 - m_Right, [343](#)
 - m_Top, [343](#)
- CCreateFilterNet, [46](#)
 - ~CCreateFilterNet, [46](#)
 - CCreateFilterNet, [46](#)
 - CutoffFrequency, [47](#)
 - FindFilter, [47](#)
 - GetBiQuad, [47](#)
 - GetBiQuads, [47](#)
 - NumCoefSets, [47](#)
 - Order, [47](#)
 - SampleRate, [47](#)
 - Scale, [47](#)
- CDacCalibrationFunctionNet, [48](#)
 - !CDacCalibrationFunctionNet, [48](#)
 - ~CDacCalibrationFunctionNet, [48](#)
 - BurnDacOffset, [49](#)
 - CDacCalibrationFunctionNet, [48](#)
 - GetDacOffset, [49](#)
 - SetDacOffset, [49](#)
- CDacqGroupChannelGenericSelectionNet, [49](#)
 - CDacqGroupChannelGenericSelectionNet, [50](#)
- CDacqGroupChannelSelectionNet, [50](#)
 - CDacqGroupChannelSelectionNet, [50](#)
- CDacqGroupChannelSelectionTemplateNet
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [51](#)
- CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [51](#)
- CDacqGroupChannelSelectionTemplateNet, [51](#)

- EnableChannelsInGroup, [51](#), [52](#)
- GetDeviceGroupChannelInfos, [52](#)
- GetEnabledChannelsInGroup, [52](#)
- GetGroupID, [52](#)
- GetGroupNumberOfChannels, [52](#)
- GetGroupSampleSize, [53](#)
- GetGroupType, [53](#)
- GetNumberOfSupportedGroups, [53](#)
- CDeviceGroupChannelInfoGenericNet, [53](#)
 - CDeviceGroupChannelInfoGenericNet, [54](#)
- CDeviceGroupChannelInfoNet, [54](#)
 - CDeviceGroupChannelInfoNet, [54](#)
- CDeviceGroupChannelInfoSCUNet, [54](#)
 - CDeviceGroupChannelInfoSCUNet, [55](#)
- CDeviceGroupChannelInfoTemplateNet
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [55](#)
- CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [55](#)
 - CDeviceGroupChannelInfoTemplateNet, [55](#)
 - GroupID, [55](#)
 - GroupType, [55](#)
 - NumberOfChannels, [56](#)
- CDeviceGroupChannelInfoW2100Net, [56](#)
 - CDeviceGroupChannelInfoW2100Net, [56](#)
- CDigOutStimulatorFunctionNet, [56](#)
 - !CDigOutStimulatorFunctionNet, [58](#)
 - ~CDigOutStimulatorFunctionNet, [57](#)
 - CDigOutStimulatorFunctionNet, [57](#)
 - ClearChannel, [58](#)
 - GetGlobalRepeat, [58](#)
 - GetNumberOfChannels, [58](#)
 - GetStartTriggerSlope, [58](#)
 - GetStopTriggerSlope, [59](#)
 - PrepareChannelData, [59](#)
 - SendChannelData, [59](#)
 - SetGlobalRepeat, [60](#)
 - SetStartTriggerSlope, [60](#)
 - SetStopTriggerSlope, [60](#)
- CEncapsulatorDeviceNet, [60](#)
 - CEncapsulatorDeviceNet, [61](#)
 - GetRoboFluidDevice, [61](#)
- CExternDTesterDeviceNet, [61](#)
 - !CExternDTesterDeviceNet, [62](#)
 - ~CExternDTesterDeviceNet, [62](#)
 - CExternDTesterDeviceNet, [62](#)
 - Read, [62](#)
 - Read2, [63](#)
 - Write, [63](#)
 - Write2, [63](#)
- CFilterCoefficientsNet, [63](#)
 - ~CFilterCoefficientsNet, [65](#)
 - A, [65](#)
 - A1, [65](#)
 - A2, [66](#)
 - B, [66](#)
 - B0, [66](#)
 - B1, [66](#)
 - B2, [66](#)
- CFilterCoefficientsNet, [64](#), [65](#)
- Equals, [65](#)
- SetAFormat, [65](#)
- SetBFormat, [65](#)
- UIntA1, [66](#)
- UIntA2, [66](#)
- UIntB0, [66](#)
- UIntB1, [66](#)
- UIntB2, [66](#)
- CFilterConfigurationNet, [67](#)
 - CFilterConfigurationNet, [67](#)
 - EraseFilterParameterPermanent, [67](#)
 - GetHighpassFilterEnable, [67](#)
 - ResetHighpassFilter, [67](#)
 - SetFilterParameter, [68](#)
 - SetFilterParameterPermanent, [68](#)
 - SetHighpassFilterEnable, [68](#)
- CFilterConfigurationRegisterNet, [68](#)
 - CFilterConfigurationRegisterNet, [69](#)
 - EraseFilterParameterPermanent, [69](#)
 - SetFilterParameter, [69](#)
 - SetFilterParameterPermanent, [69](#), [70](#)
- CFilterPropertyNet, [70](#)
 - ~CFilterPropertyNet, [70](#)
 - CFilterPropertyNet, [70](#)
 - CornerFrequency, [71](#)
 - CornerFrequencycymHz, [71](#)
 - FilterActive, [71](#)
 - FilterBand, [71](#)
 - FilterFamily, [71](#)
 - FilterType, [71](#)
 - Order, [71](#)
 - ToString, [71](#)
- CFluidControlDeviceNet, [72](#)
 - ~CFluidControlDeviceNet, [73](#)
 - CalibrateThermocouple, [73](#)
 - CFluidControlDeviceNet, [73](#)
 - GetAdc, [73](#)
 - GetDigin, [74](#)
 - GetDigout, [74](#)
 - GetPWM, [74](#)
 - GetReferenceTemperature, [74](#)
 - GetSingleValve, [75](#)
 - GetThermocoupleCalibration, [75](#)
 - GetThermocoupleNanovoltPerKelvin, [75](#)
 - GetThermocoupleTemperature, [76](#)
 - GetValve, [76](#)
 - McsBus_VoltageMode, [78](#)
 - SetDigout, [76](#)
 - SetPWM, [76](#)
 - SetSingleValve, [77](#)
 - SetThermocoupleNanovoltPerKelvin, [77](#)
 - SetValve, [77](#)
- CFYIDeviceNet, [78](#)
 - CFYIDeviceNet, [78](#)
 - FYIProgram, [79](#)
 - FYITemp, [79](#)

- Sensor, [79](#)
- CGenericDevelopDeviceNet, [79](#)
 - ~CGenericDevelopDeviceNet, [86](#)
 - CGenericDevelopDeviceNet, [86](#)
 - ClosePipe, [86](#)
 - GetBuffer, [86](#)
 - GetByteBuffer, [87](#)
 - GetIntBuffer, [87](#)
 - GetShortBuffer, [88](#)
 - GetUByteBuffer, [89](#)
 - GetUIntBuffer, [89](#)
 - GetUShortBuffer, [90](#)
 - OpenPipe, [91](#)
 - ReadPipe, [91](#)
 - ResetPipe, [91](#)
 - SetBuffer, [92](#)
 - SetByteBuffer, [92](#)
 - SetIntBuffer, [92](#)
 - SetShortBuffer, [93](#)
 - SetUByteBuffer, [94](#)
 - SetUIntBuffer, [94](#)
 - SetUShortBuffer, [95](#)
 - SetValue, [95](#)
 - WritePipe, [96](#)
- CGilsonDeviceNet, [96](#)
 - ~CGilsonDeviceNet, [97](#)
 - CGilsonDeviceNet, [97](#)
 - ConnectSlave, [97](#)
 - GetLastAnswer, [98](#)
 - m_pGilsonDevice, [98](#)
 - SendBuffered, [98](#)
 - SendImmediate, [98](#)
 - SendImmediateGetResponse, [98](#)
- ChangeSerialNumber
 - CMcsUsbFactoryNet, [209](#)
- channel
 - CMeaAudioFunctionNet::s_setaudionet, [590](#)
 - CW2100_FunctionNet::AudioChannelsNet, [27](#)
- ChannelBlock_AvailFrames
 - CMcsUsbDacqNet, [166](#)
- ChannelBlock_ReadAsFrameArrayI16
 - CMcsUsbDacqNet, [166](#), [167](#)
- ChannelBlock_ReadAsFrameArrayI32
 - CMcsUsbDacqNet, [168](#)
- ChannelBlock_ReadAsFrameArrayUI16
 - CMcsUsbDacqNet, [169](#), [170](#)
- ChannelBlock_ReadAsFrameArrayUI32
 - CMcsUsbDacqNet, [170](#), [171](#)
- ChannelBlock_ReadFramesDictI16
 - CMcsUsbDacqNet, [172](#)
- ChannelBlock_ReadFramesDictI32
 - CMcsUsbDacqNet, [172](#)
- ChannelBlock_ReadFramesDictUI16
 - CMcsUsbDacqNet, [173](#)
- ChannelBlock_ReadFramesDictUI32
 - CMcsUsbDacqNet, [173](#)
- ChannelBlock_ReadFramesI16
 - CMcsUsbDacqNet, [174](#), [175](#)
- ChannelBlock_ReadFramesI32
 - CMcsUsbDacqNet, [175](#), [176](#)
- ChannelBlock_ReadFramesUI16
 - CMcsUsbDacqNet, [177](#)
- ChannelBlock_ReadFramesUI32
 - CMcsUsbDacqNet, [178](#), [179](#)
- ChannelDataEvent
 - CMcsUsbDacqNet, [205](#)
- ChannelReset
 - CMultiBatteryChargerDeviceNet, [283](#)
- Charge
 - BatteryState, [28](#)
- ChargeRegionString
 - BatteryState, [28](#)
- ChargeString
 - BatteryState, [28](#)
- CHiClampDeviceNet, [98](#)
 - CHiClampDeviceNet, [99](#)
 - RoboDacq, [99](#)
- CHLADacqNet, [99](#)
 - CHLADacqNet, [100](#)
- CHLADeviceNet, [100](#)
 - CHLADeviceNet, [100](#)
 - HLADacq, [101](#)
 - SerialPort, [101](#)
- CHWInfo
 - CMcsUsbDacqNet::CHWInfo, [101](#)
- CIntanMea_FunctionNet, [103](#)
 - AmplifierSettle, [104](#)
 - BeginImpedanceCheck, [104](#)
 - CIntanMea_FunctionNet, [104](#)
 - GetDSPHighPassByIndex, [104](#)
 - GetImpedanceResult, [104](#)
 - GetIntanRegister, [104](#)
 - GetLowerFrequencyByIndex, [104](#)
 - GetUpperFrequencyByIndex, [105](#)
 - SetBandwidthByIndex, [105](#)
 - SetDiagnosticMode, [105](#)
 - SetDSPHighPassByIndex, [105](#)
 - SetIntanRegister, [105](#)
- CInterfaceboardFunctionNet, [105](#)
 - !CInterfaceboardFunctionNet, [106](#)
 - ~CInterfaceboardFunctionNet, [106](#)
 - CInterfaceboardFunctionNet, [106](#)
 - GetCardinalDacqSamplerate, [107](#)
 - GetCardinalStgOutputrate, [107](#)
 - SetCardinalDacqSamplerate, [107](#)
 - SetCardinalStgOutputrate, [107](#)
- ClampAmpRestart
 - CRoboDacqNet, [357](#)
- ClearBuffers
 - CMcsUsbDacqNet, [179](#)
- ClearChannel
 - CDigOutStimulatorFunctionNet, [58](#)
- ClearChannel_PrepAndSendData
 - CStg200xDownloadNet, [455](#)
 - CStimulusFunctionNet, [462](#)
- ClearChannelData

- CStg200xDownloadBasicNet, 446
- CStimulusFunctionNet, 462
- CW2100_StimulatorFunctionNet, 508
- ClearMultiplexedData
 - CStimulusFunctionNet, 462
- ClearSTGOutput
 - CCMOSMea_FunctionNet, 34
- ClearStimulusParametersCache
 - CW2100_FunctionNet, 501
- ClearSyncData
 - CStg200xDownloadBasicNet, 446
 - CStimulusFunctionNet, 462
- ClearTable
 - CRetinalLedDeviceNet, 345
- ClearTableName
 - CWarnerValveControllerDeviceNet, 536
- ClearUserDefinedNameCache
 - CW2100_FunctionNet, 501
- ClearValveTable
 - CWarnerValveControllerDeviceNet, 536
- CLIH3DeviceNet, 107
 - !CLIH3DeviceNet, 109
 - ~CLIH3DeviceNet, 109
 - CLIH3DeviceNet, 109
 - DummyCommand, 109
 - EnableUserTrigger, 110
 - GetADCOffset, 110
 - GetDaclIdleValue, 110
 - GetDacqRunStatus, 110
 - GetDacUseIdleValue, 111
 - GetDigInState, 111
 - GetEEpromPage, 111
 - GetSampleInterval, 112
 - IsUserTriggerEnabled, 112
 - ReadClipping, 112
 - SendCommand, 112
 - SetADCOffset, 113
 - SetDaclIdleValue, 113
 - SetDacUseIdleValue, 113
 - SetDigOutState, 113
 - SetEEpromPage, 114
 - SetSampleInterval, 114
 - StimulusFunction, 114
- CloseAllValves
 - CRoboFluidDeviceNet, 382
- ClosePipe
 - CGenericDevelopDeviceNet, 86
- ClosePlateClamp
 - CMultiwellDeviceNet, 293
- CMcsBus_AxisParametersNet, 114
 - ~CMcsBus_AxisParametersNet, 115
 - CMcsBus_AxisParametersNet, 115
 - GetAxisParametersSignedEeprom, 115
 - GetAxisParametersUnsignedEeprom, 115
 - SetAxisParametersEeprom, 115, 116
- CMcsBus_ExtensionNet, 116
 - ~CMcsBus_ExtensionNet, 116
 - CMcsBus_ExtensionNet, 116
- GetLEDSwitch, 117
- SetLEDSwitch, 117
- CMcsBus_FYIExtensionNet, 117
 - ~CMcsBus_FYIExtensionNet, 117
 - CMcsBus_FYIExtensionNet, 117
- GetDIO, 118
- GetSingleHeater, 118
- GetValves, 118
- SetDIO, 118
- SetSingleHeater, 118
- SetValves, 118
- CMcsBus_MotorControlNet, 119
 - ~CMcsBus_MotorControlNet, 122
 - CMcsBus_MotorControlNet, 122
 - GetMCAcceleration, 122
 - GetMCAccelerationEeprom, 122
 - GetMCAccelerationShortCommand, 122
 - GetMCAXisRevisionEeprom, 122
 - GetMCBreakCurrent, 122
 - GetMCBreakCurrentEeprom, 123
 - GetMCConfig, 123
 - GetMCConfigEeprom, 123
 - GetMCCurrent, 123
 - GetMCCurrentEeprom, 123
 - GetMCCurrentMode, 123
 - GetMCCurrentModeEeprom, 123
 - GetMCCurrentModeShortCommand, 124
 - GetMCCurrentPosition, 124
 - GetMCCurrentShortCommand, 124
 - GetMCCurrentSpeed, 124
 - GetMCMaxAcceleration, 124
 - GetMCMaxAccelerationEeprom, 124
 - GetMCMaxCurrent, 124
 - GetMCMaxCurrentEeprom, 125
 - GetMCMaxSpeed, 125
 - GetMCMaxSpeedEeprom, 125
 - GetMCMaxTravel, 125
 - GetMCMaxTravelEeprom, 125
 - GetMCMaxTravelShortCommand, 125
 - GetMCMovement, 125
 - GetMCNewPosition, 126
 - GetMCOutputOnOff, 126
 - GetMCPhase, 126
 - GetMCPhaseOffset, 126
 - GetMCReference, 126
 - GetMCReferenceCurrent, 126
 - GetMCReferenceCurrentEeprom, 126
 - GetMCRegulatorGain, 127
 - GetMCRegulatorGainEeprom, 127
 - GetMCScalingFactor, 127
 - GetMCScalingFactorEeprom, 127
 - GetMCSpeed, 127
 - GetMCSpeedEeprom, 127
 - GetMCSpeedShortCommand, 127
 - GetMCSpeedUnitEeprom, 128
 - GetMCStandbyCurrent, 128
 - GetMCStandbyCurrentEeprom, 128
 - GetMCStandbyTime, 128

- GetMCStandbyTimeEeprom, 128
- GetSubChannel, 128
- SetMCAcceleration, 128
- SetMCAccelerationEeprom, 129
- SetMCAccelerationShortCommand, 129
- SetMCAxisRevisionEeprom, 129
- SetMCBreakCurrent, 129
- SetMCBreakCurrentEeprom, 129
- SetMCConfig, 129
- SetMCConfigEeprom, 130
- SetMCCurrent, 130
- SetMCCurrentEeprom, 130
- SetMCCurrentMode, 130
- SetMCCurrentModeEeprom, 130
- SetMCCurrentModeShortCommand, 130
- SetMCCurrentPosition, 131
- SetMCCurrentShortCommand, 131
- SetMCMaxAcceleration, 131
- SetMCMaxAccelerationEeprom, 131
- SetMCMaxCurrent, 131
- SetMCMaxCurrentEeprom, 131
- SetMCMaxSpeed, 132
- SetMCMaxSpeedEeprom, 132
- SetMCMaxTravel, 132
- SetMCMaxTravelEeprom, 132
- SetMCMaxTravelShortCommand, 132
- SetMCNewPosition, 132
- SetMCOutputOnOff, 133
- SetMCReference, 133
- SetMCReferenceCurrent, 133
- SetMCReferenceCurrentEeprom, 133
- SetMCRegulatorGain, 133
- SetMCRegulatorGainEeprom, 133
- SetMCRotation, 134
- SetMCScalingFactor, 134
- SetMCScalingFactorEeprom, 134
- SetMCSpeed, 134
- SetMCSpeedEeprom, 134
- SetMCSpeedShortCommand, 134
- SetMCSpeedUnitEeprom, 135
- SetMCStandbyCurrent, 135
- SetMCStandbyCurrentEeprom, 135
- SetMCStandbyTime, 135
- SetMCStandbyTimeEeprom, 135
- SetSubChannel, 135
- StartMCMovement, 136
- StopMCMovement, 136
- CMcsBus_SensorNet, 136
 - ~CMcsBus_SensorNet, 138
 - CatchAmpGetAdcMean, 138
 - CatchAmpGetAdcValue, 138
 - CatchAmpGetAdcValueH, 138
 - CatchAmpGetAdcValueL, 138
 - CatchAmpGetDacAmplitude, 139
 - CatchAmpGetDacEnable, 139
 - CatchAmpGetDacOffset, 139
 - CatchAmpGetPwmEnable, 139
 - CatchAmpSetDacAmplitude, 139
 - CatchAmpSetDacEnable, 139
 - CatchAmpSetDacOffset, 139
 - CatchAmpSetPwmEnable, 139
 - CMcsBus_SensorNet, 138
 - Get2AnalogInput, 140
 - Get2DigitalInput, 140
 - Get4ADC, 140
 - Get4ADCAverage, 140
 - Get4ADCCatchampAverageShift, 140
 - Get4ADCMode, 140
 - Get4DAC, 140
 - GetADCs, 140
 - GetADCsLoop, 141
 - GetBubbleStatus, 141
 - GetDACs, 141
 - GetDetectionThreshold, 141
 - GetDetectorValue, 141
 - GetLatency, 141
 - GetLatencyCounter, 141
 - GetMinimalThreshold, 141
 - GetMovePump, 142
 - GetPiezoState, 142
 - GetPressure, 142
 - GetPressureOffset, 142
 - GetRegulationTimeouts, 142
 - GetRegulatorFactor, 143
 - GetRegulatorOnOff, 143
 - GetRegulatorStatus, 143
 - GetRotatePump, 143
 - GetSamplePeriode, 143
 - GetSollPressure, 143
 - GetSyncState, 143
 - Set4ADCCatchampAverageShift, 144
 - Set4ADCMode, 144
 - Set4DAC, 144
 - SetDACs, 144
 - SetDetectionThreshold, 144
 - SetLatency, 144
 - SetMinimalThreshold, 144
 - SetMovePump, 145
 - SetPiezoState, 145
 - SetPressureOffset, 145
 - SetRegulationTimeouts, 145
 - SetRegulatorFactor, 145
 - SetRegulatorOnOff, 145
 - SetRotatePump, 145
 - SetSamplePeriode, 146
 - SetSollPressure, 146
 - StartSync, 146
 - TactSwitchGetState, 146
 - TactSwitchSetDisplay, 146
- CMcsBus_TempSensorNet, 146
 - ~CMcsBus_TempSensorNet, 147
 - CMcsBus_TempSensorNet, 147
 - GetNanoVoltsPerKelvin, 147
 - GetTemperatur, 147
 - GetThermoOffset, 147
 - GetThermoTemp, 148

- GetThermoVoltage, 148
- SetNanoVoltsPerKelvin, 148
- SetThermoOffset, 148
- CMcsBus_VoltageModeNet, 148
 - ~CMcsBus_VoltageModeNet, 149
 - CMcsBus_VoltageModeNet, 149
 - GetVMMMaxNegativeCurrent, 150
 - GetVMMMaxNegativeCurrentEeprom, 150
 - GetVMMMaxNegativeVoltage, 150
 - GetVMMMaxNegativeVoltageEeprom, 150
 - GetVMMMaxPositiveCurrent, 150
 - GetVMMMaxPositiveCurrentEeprom, 150
 - GetVMMMaxPositiveVoltage, 150
 - GetVMMMaxPositiveVoltageEeprom, 151
 - GetVMOOutputOnOff, 151
 - GetVMVoltage, 151
 - SetVMMMaxNegativeCurrent, 151
 - SetVMMMaxNegativeCurrentEeprom, 151
 - SetVMMMaxNegativeVoltage, 151
 - SetVMMMaxNegativeVoltageEeprom, 151
 - SetVMMMaxPositiveCurrent, 152
 - SetVMMMaxPositiveCurrentEeprom, 152
 - SetVMMMaxPositiveVoltage, 152
 - SetVMMMaxPositiveVoltageEeprom, 152
 - SetVMOOutputOnOff, 152
 - SetVMVoltage, 152
- CMcsBusNet, 153
 - ~CMcsBusNet, 153
 - CMcsBusNet, 153
 - CMcsBusNet::GetMode, 154
 - CMcsBusNet::GetModeEeprom, 154
 - CMcsBusNet::SetMode, 154
 - CMcsBusNet::SetModeEeprom, 154
 - GetBusAddress, 154
 - GetBusAddressEeprom, 154
 - GetCommand, 154, 155
 - GetHWRevisionEeprom, 155
 - SetBusAddress, 155
 - SetBusAddressEeprom, 155
 - SetCommand, 156
 - SetHWRevisionEeprom, 156
- CMcsBusNet::GetMode
 - CMcsBusNet, 154
- CMcsBusNet::GetModeEeprom
 - CMcsBusNet, 154
- CMcsBusNet::SetMode
 - CMcsBusNet, 154
- CMcsBusNet::SetModeEeprom
 - CMcsBusNet, 154
- CMcsUsbDacqNet, 157
 - ~CMcsUsbDacqNet, 163
 - AddSelectedChannelsQueue, 163–165
 - ChannelBlock_AvailFrames, 166
 - ChannelBlock_ReadAsFrameArrayI16, 166, 167
 - ChannelBlock_ReadAsFrameArrayI32, 168
 - ChannelBlock_ReadAsFrameArrayUI16, 169, 170
 - ChannelBlock_ReadAsFrameArrayUI32, 170, 171
 - ChannelBlock_ReadFramesDictI16, 172
 - ChannelBlock_ReadFramesDictI32, 172
 - ChannelBlock_ReadFramesDictUI16, 173
 - ChannelBlock_ReadFramesDictUI32, 173
 - ChannelBlock_ReadFramesI16, 174, 175
 - ChannelBlock_ReadFramesI32, 175, 176
 - ChannelBlock_ReadFramesUI16, 177
 - ChannelBlock_ReadFramesUI32, 178, 179
 - ChannelDataEvent, 205
 - ClearBuffers, 179
 - CMcsUsbDacqNet, 163
 - CMcsUsbDacqNet::GetFilterProperties, 179
 - Error_Callback_Aquisition_Stopped, 204
 - Error_Callback_Data_lost, 204
 - Error_Callback_Frames_Lost, 204
 - Error_Callback_Packet_Error, 204
 - Error_Callback_Queue_Full, 204
 - Error_Callback_RingQueue_Full, 204
 - ErrorEvent, 205
 - GetAdapterType, 180
 - GetAdcDataFormat, 180
 - GetAdcZero, 180
 - GetAnalogValueUnit, 180
 - GetChannelDataFillSize, 180
 - GetChannelLayout, 180
 - GetChannelsInBlock, 181
 - GetDataFormat, 181
 - GetDataMode, 181
 - GetDigitalSource, 181–183
 - GetFilterProperty, 183
 - GetGroupChannelDataI16, 183
 - GetGroupChannelDataI32, 184
 - GetGroupChannelDataUI16, 184
 - GetGroupChannelDataUI32, 185
 - GetHardwareMaxRange, 186
 - GetHardwareMinRange, 186
 - GetMaxSamplingFrequency, 186
 - GetMeaLayout, 186
 - GetMinSamplingFrequencyStepsize, 186
 - GetNumberOfDataBits, 187
 - GetPoti, 187
 - GetResolutionPerDigit, 187
 - GetSamplerate, 187
 - GetVoltageRangeIndex, 187
 - GetVoltageRangeInMicroVolt, 187
 - GetVoltageRangeInMilliVolt, 188
 - HWInfo, 188
 - Samplerate, 204
 - SendStartDacq, 188
 - SendStartStgAndDacq, 188
 - SendStopDacq, 189
 - SendStopStgAndDacq, 189
 - SendStopStgAndDacqWithOptions, 189
 - SetDataMode, 190
 - SetDigitalSource, 190–192
 - SetPoti, 192
 - SetSamplerate, 192
 - SetSelectedChannels, 192–194
 - SetSelectedChannelsQueue, 195–197

- SetSelectedData, 197–199
- SetupGroupDacqQueue, 199
- SetVoltageRangeByIndex, 199
- SetVoltageRangeInMicroVolt, 199
- StartDacq, 199–201
- StartLoop, 201, 202
- StopDacq, 203
- StopLoop, 204
- CMcsUsbDacqNet::CHWInfo, 101
 - CHWInfo, 101
 - GetAvailableSampleRates, 102
 - GetAvailableVoltageRangesInMicroVolt, 102
 - GetAvailableVoltageRangesInMicroVoltAnd-StringsInMilliVolt, 102
 - GetNumberOfHWADCCChannels, 102
 - GetNumberOfHWDigitalChannels, 102
 - IsDigitalChannelDedicated, 103
- CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, 499
 - CVoltageRangeInfoNet, 499
 - VoltageRangeDisplayStringMilliVolt, 499
 - VoltageRangeInMicroVolt, 499
- CMcsUsbDacqNet::GetFilterProperties
 - CMcsUsbDacqNet, 179
- CMcsUsbDeviceStatePushFunctionNet, 205
 - CMcsUsbDeviceStatePushFunctionNet, 205
 - McsUsbDeviceStateEvent, 206
 - TriggerStatus, 206
- CMcsUsbDeviceStatePushNet, 206
 - CMcsUsbDeviceStatePushNet, 206
 - McsUsbDeviceStateEvent, 207
 - TriggerStatus, 207
- CMcsUsbFactoryNet, 207
 - ~CMcsUsbFactoryNet, 209
 - ChangeSerialNumber, 209
 - CMcsUsbFactoryNet, 209
 - Coldstart, 209
 - CompareFirmware, 209
 - DownloadFirmware, 209
 - FindFirmwareVersionMagicInBuffer, 209
 - FX3MCSDDataAddress, 215
 - FX3MCSDDataDeviceIdOffset, 215
 - FX3MCSDDataVersionOffset, 215
 - GetChecksumFromFX3Image, 210
 - GetDestination, 210
 - GetDestinationDisplayLabel, 210
 - GetDestinationName, 210
 - GetDestinationSerialNumber, 210
 - GetDestinationTargetAddress, 210
 - GetFirmwareVersionFromFile, 211
 - GetFirmwareVersionFromHexFile, 211
 - GetNumDestinations, 211
 - GetUSBDeviceIDFromFX3Image, 211
 - GetUserCodeFromBitFile, 211
 - GetUserCodeFromFlash, 211
 - GetXilinxFlashOffset, 212
 - GetXilinxFlashReadCommand, 212
 - LoadUserFirmware, 212
 - ReadBlockFromFlash, 212
 - ReadBlockFromIFBGlobalEEPROM, 213
 - ReadBlockFromNVMEM, 213
 - SetDestinationSerialNumber, 213
 - UpdateFirmware, 213, 214
- CMcsUsbFunctionNet, 215
 - !CMcsUsbFunctionNet, 216
 - ~CMcsUsbFunctionNet, 216
 - CMcsUsbFunctionNet, 216
 - m_pMcsUsb, 217
 - m_pMcsUsbFunction, 217
 - ThrowCUsbExceptionNetOnError, 216
- CMcsUsbFunctionPointerContainer, 217
- CMcsUsbListEntryNet, 217
 - ~CMcsUsbListEntryNet, 218
 - DeviceId, 220
 - DeviceName, 220
 - Equals, 218
 - GetEntry, 218, 219
 - GetEntryCount, 219
 - HwVersion, 220
 - Manufacturer, 220
 - Product, 221
 - SerialNumber, 221
 - SetStringFormat, 220
 - ToString, 220
- CMcsUsbListNet, 221
 - !CMcsUsbListNet, 222
 - ~CMcsUsbListNet, 222
 - CMcsUsbListNet, 222
 - Count, 223
 - DeviceArrival, 223
 - DeviceRemoval, 223
 - GetNumberOfDevices, 222
 - GetUsbListEntries, 222
 - GetUsbListEntry, 222
 - IsDeviceTypeOf, 223
 - SetStringFormat, 223
- CMcsUsbNet, 224
 - !CMcsUsbNet, 228
 - ~CMcsUsbNet, 228
 - AddSoftwareKey, 228
 - AssociateToThis, 228
 - CMcsUsbNet, 228
 - Connect, 228, 229
 - Disconnect, 230
 - EmptyKey, 230
 - EnableExceptions, 230
 - EraseEEPROMRegisterPreconfig, 230
 - GetConfiguration, 231
 - GetDeviceCannotStallOutRequests, 231
 - GetDeviceCapableSpeed, 231
 - GetDeviceEnum, 231
 - GetDeviceId, 231
 - GetDeviceRootHubVendorEnum, 231
 - GetDeviceRootHubVendorID, 231
 - GetDeviceRootHubVendorName, 231
 - GetDeviceSpeed, 232

- GetErrorText, [232](#)
- GetFirmwareVersion, [232](#)
- GetHardwareRevision, [232](#)
- GetHeadstageActive, [233](#)
- GetHeadstageID, [233](#)
- GetHeadstagePresent, [233](#)
- GetIdent, [233](#)
- GetLastUSBError, [234](#)
- GetMea21UsbPort, [234](#)
- GetNumConfigurations, [234](#)
- GetSerialNumber, [234](#)
- GetSoftwareKey, [234](#)
- GetSoftwareKeyString, [234](#)
- GetStatus, [234](#)
- GetStatusOfLastCommand, [235](#)
- GetUsbListEntry, [235](#)
- GetVersion, [235](#)
- HasSoftwareKey, [235](#)
- IsConnected, [235](#)
- IsDeviceHighSpeed, [235](#)
- IsDeviceHighSpeedCapable, [236](#)
- IsExceptionsEnabled, [236](#)
- MultibootGetCypressImageId, [236](#)
- MultibootGetImageId, [236](#)
- MultibootGetSelectedImage, [236](#)
- MultibootSelectImage, [236](#)
- ReadEepromRegisterPreconfig, [237](#)
- ReadRegister, [237](#)
- ReadRegister32, [237](#)
- ReadRegisterTimeSlot, [237](#)
- RemoveSoftwareKey, [238](#)
- RescanHeadstage, [238](#)
- SerialNumber, [245](#)
- SetConfiguration, [238](#)
- SetSoftwareKey, [238](#)
- Status_AlreadyConfigured, [240](#)
- Status_BadStartFrame, [241](#)
- Status_Btstuff, [241](#)
- Status_BufferOverrun, [241](#)
- Status_BufferUnderrun, [241](#)
- Status_Canceled, [241](#)
- Status_Canceling, [241](#)
- Status_ConnectedPipes, [241](#)
- Status_ControlNotOwned, [241](#)
- Status_Crc, [241](#)
- Status_DataOverrun, [241](#)
- Status_DataToggleMismatch, [241](#)
- Status_DataUnderrun, [242](#)
- Status_DeviceLocked, [242](#)
- Status_DeviceNotFound, [242](#)
- Status_DeviceRemoved, [242](#)
- Status_DevNotResponding, [242](#)
- Status_EndpointHalted, [242](#)
- Status_ErrorBusy, [242](#)
- Status_ErrorShortTransfer, [242](#)
- Status_Fifo, [242](#)
- Status_FrameControlOwned, [242](#)
- Status_InternalHcError, [242](#)
- Status_InvalidParameter, [243](#)
- Status_InvalidPipeHandle, [243](#)
- Status_InvalidUrbFunction, [243](#)
- Status_IoPending, [243](#)
- Status_IoTimeout, [243](#)
- Status_IsochRequestFailed, [243](#)
- Status_LastUsbErrorMismatch, [243](#)
- Status_NoBandwidth, [243](#)
- Status_NoMemory, [243](#)
- Status_NoSuchDevice, [243](#)
- Status_NotAccessed, [243](#)
- Status_NotSupported, [244](#)
- Status_PidCheckFailure, [244](#)
- Status_PipeNotLinked, [244](#)
- Status_RequestFailed, [244](#)
- Status_RequestMutexFailed, [244](#)
- Status_RequestMutexTimeout, [244](#)
- Status_Stall, [244](#)
- Status_Unconfigured, [244](#)
- Status_UnexpectedPid, [244](#)
- ThrowCUsbExceptionNetOnError, [238](#)
- TxnGetSerialNumber, [238](#)
- TxnSetSerialNumber, [238](#)
- TxnTestMemoryReadAndCheck, [238](#)
- TxnTestMemoryWrite, [239](#)
- ValidKey, [239](#)
- WPAError_ScanningIsPending, [244](#)
- WriteEepromRegisterPreconfig, [239](#)
- WriteRegister, [239, 240](#)
- WriteRegister32, [240](#)
- WriteRegisterArray, [240](#)
- WriteRegisterTimeSlot, [240](#)
- WriteRegisterValue, [240](#)
- CMcsUsbPointerContainer, [245](#)
- CMEA2100x256FunctionNet, [245](#)
 - !CMEA2100x256FunctionNet, [246](#)
 - ~CMEA2100x256FunctionNet, [246](#)
 - CMEA2100x256FunctionNet, [245, 246](#)
 - GetLayoutConfiguration, [246](#)
 - SetLayoutConfiguration, [246](#)
- CMeaAudioFunctionNet, [247](#)
 - CMeaAudioFunctionNet, [247](#)
 - GetAudioChannels, [248, 249](#)
 - GetNumberOfAudioChannels, [249](#)
 - SetAudioChannels, [249, 250](#)
- CMeaAudioFunctionNet::s_setaudionet, [590](#)
 - amplification, [590](#)
 - channel, [590](#)
- CMeaCleanDeviceNet, [250](#)
 - !CMeaCleanDeviceNet, [252](#)
 - ~CMeaCleanDeviceNet, [251](#)
 - CMeaCleanDeviceNet, [251](#)
 - GetCycle, [252](#)
 - GetCycles, [252](#)
 - GetMaxVoltage, [252](#)
 - GetMinVoltage, [252](#)
 - GetOutputVoltage, [252](#)
 - GetSlope, [253](#)

- IsRunning, [253](#)
- SetCycles, [253](#)
- SetMaxVoltage, [253](#)
- SetMinVoltage, [254](#)
- SetSlope, [254](#)
- Start, [254](#)
- Stop, [254](#)
- CMeaCoatDeviceNet, [254](#)
 - !CMeaCoatDeviceNet, [256](#)
 - ~CMeaCoatDeviceNet, [256](#)
 - CMeaCoatDeviceNet, [256](#)
 - GetCurrentCycle, [256](#)
 - GetCycles, [256](#)
 - GetDuration, [256](#)
 - GetMaxCurrent, [256](#)
 - GetOffsetCurrent, [257](#)
 - GetOutputCurrent, [257](#)
 - GetPauseDuration, [257](#)
 - GetSlope, [257](#)
 - GetTimeInPause, [257](#)
 - GetTimeInPlateau, [258](#)
 - IsRunning, [258](#)
 - SetCycles, [258](#)
 - SetDuration, [258](#)
 - SetMaxCurrent, [259](#)
 - SetOffsetCurrent, [259](#)
 - SetPauseDuration, [259](#)
 - SetSlope, [259](#)
 - Start, [259](#)
 - Stop, [260](#)
- CMeaDeviceNet, [260](#)
 - ~CMeaDeviceNet, [262](#)
 - AnalogGain, [267](#)
 - CMeaDeviceNet, [261](#)
 - EnableChecksum, [262](#)
 - EnableDigitalIn, [262](#), [263](#)
 - EnableTimestamp, [263](#)
 - Gain, [267](#)
 - GetAnalogGain, [264](#)
 - GetEnumerationSpeed, [264](#)
 - GetGain, [264](#)
 - MeaAudioFunctionNet, [267](#)
 - MeaDigitalDataFunctionNet, [267](#)
 - MeaFeedbackFunctionNet, [267](#)
 - MeFunctionNet, [268](#)
 - SetDigitalOut, [264](#)
 - SetNumberOfAnalogChannels, [264](#)
 - SetNumberOfChannels, [265](#), [266](#)
 - SetTriggerMaskValue, [266](#)
 - SetTriggerPeriod, [267](#)
 - W2100_FunctionNet, [268](#)
 - WClassicFunctionNet, [268](#)
- CMeaDigitalDataFunctionNet, [268](#)
 - CMeaDigitalDataFunctionNet, [268](#)
 - GetDigitalData, [269](#)
 - SetDigitalData, [269](#)
- CMeaFeedbackFunctionNet, [270](#)
 - CMeaFeedbackFunctionNet, [271](#)
 - FeedbackGetSampleTimerCount, [271](#)
 - FeedbackSetAnalogSource, [271](#)
 - FeedbackSetChannelFilter, [271](#)
 - FeedbackSetDigitalMapping, [271](#)
 - FeedbackSetFeedback, [271](#)
 - FeedbackSetFilterOff, [272](#)
 - FeedbackSetFilterParameter, [272](#)
 - FeedbackSetFilterParameter32, [272](#)
 - FeedbackSetGlobalChannelFilter, [272](#)
 - FeedbackSetIIRFilterParameter, [272](#)
 - FeedbackSetLogic, [272](#)
 - FeedbackSetMkFilter, [272](#)
 - FeedbackSetNumberOfLogics, [273](#)
 - FeedbackSetNumberOfRateCounter, [273](#)
 - FeedbackSetNumberOfRateDetectors, [273](#)
 - FeedbackSetNumberOfSpikeDetectors, [273](#)
 - FeedbackSetNumberOfTriggers, [273](#)
 - FeedbackSetRateCounter, [273](#)
 - FeedbackSetRateDetector, [273](#)
 - FeedbackSetSpikeDetectorThreshold, [273](#)
 - FeedbackSetTrigger, [274](#)
- CMealImpedanceDeviceNet, [274](#)
 - ~CMealImpedanceDeviceNet, [274](#)
 - CMealImpedanceDeviceNet, [274](#)
 - GetAdapterCode, [275](#)
 - GetArraySize, [275](#)
 - GetImpedanceTestFrequency, [275](#)
 - GetReady, [275](#)
 - GetResult, [275](#)
 - SetImpedanceTestFrequency, [275](#)
 - StartMeasurement, [275](#)
- CMeasureTableDeviceNet, [275](#)
 - CMeasureTableDeviceNet, [276](#)
 - Sensor, [276](#)
- CMeaSwitchDeviceNet, [276](#)
 - ~CMeaSwitchDeviceNet, [277](#)
 - CMeaSwitchDeviceNet, [277](#)
 - GetNumber, [277](#)
 - GetPattern, [277](#)
 - GetPatternBool, [278](#)
 - SetPattern, [278](#)
 - SetPatternBool, [278](#)
- CMeaUSBDeviceNet, [278](#)
 - ~CMeaUSBDeviceNet, [279](#)
 - CMeaUSBDeviceNet, [279](#)
- CMeFunctionNet, [279](#)
 - !CMeFunctionNet, [280](#)
 - ~CMeFunctionNet, [280](#)
 - CMeFunctionNet, [280](#)
 - SetTransformer, [280](#)
- CMosMea
 - CCMOSMeaDeviceNet, [45](#)
- CMultiBatteryChargerDeviceNet, [281](#)
 - !CMultiBatteryChargerDeviceNet, [282](#)
 - ~CMultiBatteryChargerDeviceNet, [282](#)
 - CapacityTest, [282](#)
 - ChannelReset, [283](#)
 - CMultiBatteryChargerDeviceNet, [282](#)

- GetBatteryVoltage, [283](#)
- GetChannels, [283](#)
- GetChannelState, [283](#)
- GetChargeCapacity, [284](#)
- GetChargeCurrent, [284](#)
- GetChargingMode, [284](#)
- GetChargingPCoefficient, [284](#)
- GetDischargeCapacity, [285](#)
- GetDischargeCurrent, [285](#)
- GetDischargeCurrentSetPoint, [285](#)
- GetFinalDischargeVoltage, [286](#)
- GetRatedCapacity, [286](#)
- SetChargingMode, [286](#)
- SetChargingPCoefficient, [286](#)
- SetDischargeCurrentSetPoint, [287](#)
- SetFinalDischargeVoltage, [287](#)
- SetRatedCapacity, [287](#)
- SetRatedCapacityVolatile, [287](#)
- CMultiwellCallbackFunctionNet, [288](#)
 - !CMultiwellCallbackFunctionNet, [289](#)
 - ~CMultiwellCallbackFunctionNet, [289](#)
- CMultiwellCallbackFunctionNet, [289](#)
 - GetPlateClampStateByHeadstage, [289](#)
 - GetPlateClampStateByHeadstageEvent, [291](#)
 - GetPlateTypeByHeadstage, [290](#)
 - GetPlateTypeByHeadstageEvent, [291](#)
 - IsPlateTypeValidByHeadstage, [290](#)
 - IsPlateTypeValidByHeadstageEvent, [291](#)
 - OnGetPlateClampStateByHeadstage, [290](#)
 - OnGetPlateTypeByHeadstage, [290](#)
 - OnIsPlateTypeValidByHeadstage, [290](#)
- CMultiwellDeviceNet, [291](#)
 - !CMultiwellDeviceNet, [293](#)
 - ~CMultiwellDeviceNet, [293](#)
- ClosePlateClamp, [293](#)
- CMultiwellDeviceNet, [293](#)
 - GetPlateClampLockState, [293](#)
 - GetPlateClampState, [293](#)
 - GetPlateClampStateByHeadstage, [293](#)
 - GetPlateMux, [294](#)
 - GetPlateMuxByHeadstage, [294](#)
 - GetPlateType, [294](#)
 - GetPlateTypeByHeadstage, [294](#)
 - IsPlateTypeValid, [295](#)
 - IsPlateTypeValidByHeadstage, [295](#)
 - LockPlateClamp, [295](#)
 - OpenPlateClamp, [295](#)
 - SetPlateMux, [295](#)
 - SetPlateMuxByHeadstage, [296](#)
 - SetPlateType, [296](#)
 - SetPlateTypeByHeadstage, [296](#)
 - StopPlateClamp, [296](#)
 - UnlockPlateClamp, [296](#)
- CMultiwellOptoStimFunctionNet, [297](#)
 - !CMultiwellOptoStimFunctionNet, [298](#)
 - ~CMultiwellOptoStimFunctionNet, [298](#)
- CMultiwellOptoStimFunctionNet, [297](#), [298](#)
 - GetAbsMaxCurrentInMicroAmp, [298](#)
 - GetColorRgb, [298](#)
 - GetColorStr, [299](#)
 - GetMaxDurationHighCurrentInMicroSec, [299](#)
 - GetMaxDutyCycleHighCurrent, [299](#)
 - GetPermanentCurrentInMicroAmp, [299](#)
 - GetWaveLengthInNanometer, [300](#)
 - SetAbsMaxCurrentInMicroAmp, [300](#)
 - SetColorRgb, [300](#)
 - SetColorStr, [300](#)
 - SetMaxDurationHighCurrentInMicroSec, [301](#)
 - SetMaxDutyCycleHighCurrent, [301](#)
 - SetPermanentCurrentInMicroAmp, [301](#)
 - SetWaveLengthInNanometer, [301](#)
- CNF_GenDeviceNet, [302](#)
 - ~CNF_GenDeviceNet, [302](#)
- CNF_GenDeviceNet, [302](#)
 - Set_Values, [302](#)
- COctoPotDeviceNet, [302](#)
 - BurnAdcOffset, [303](#)
 - BurnDacOffset, [304](#)
 - COctoPotDeviceNet, [303](#)
 - EnableChecksum, [304](#)
 - EnableDigitalIn, [304](#)
 - EnableTimestamp, [304](#)
 - GetAdcOffset, [304](#)
 - GetDacOffset, [304](#)
 - PatternListStart, [304](#)
 - RampStart, [304](#)
 - ResetAdcOffset, [304](#)
 - ResetDacOffset, [304](#)
 - SetAdcOffset, [305](#)
 - SetAmplificationSwitch, [305](#)
 - SetBathclamp, [305](#)
 - SetChannelSwitch, [305](#)
 - SetDacAutoControl, [305](#)
 - SetDacOffset, [305](#)
 - SetDacValue, [305](#)
 - SetNumberOfChannels, [305](#)
 - SetOutputRate, [305](#)
 - SetPatternListEntry, [306](#)
 - SetPidParameter, [306](#)
 - SetRampParameter, [306](#)
 - SetSineParameter, [306](#)
 - SineStart, [306](#)
- CokuvisionStimulatorDeviceNet, [306](#)
 - ~CokuvisionStimulatorDeviceNet, [307](#)
- CokuvisionStimulatorDeviceNet, [307](#)
 - GetCheckVoltage, [307](#)
 - GetCurrentFactor, [307](#)
 - GetDACOffset, [308](#)
 - GetMaxPower, [308](#)
 - GetMaxVoltage, [308](#)
 - GetPulseform, [308](#)
 - GetRTC, [308](#)
 - GetStimulatorStatus, [308](#)
 - GetVoltage, [308](#)
 - SetCheckVoltage, [309](#)
 - SetCurrentFactor, [309](#)

- SetDACOffset, 309
- SetMaxPower, 309
- SetMaxVoltage, 309
- SetPulseform, 309
- SetRTC, 309
- Coldstart
 - CMcsUsbFactoryNet, 209
- CompareFirmware
 - CMcsUsbFactoryNet, 209
- CompareTo
 - HeadStageIDType, 576
- CompensateElectrodeOffset
 - CWarnerUssingFunctionNet, 518
- Connect
 - CMcsUsbNet, 228, 229
 - CRFFunctionNet, 347
- ConnectDevice
 - CRadioControlledDevicesNet, 342
- ConnectedImp
 - CPositionImpDeviceNet, 321
- ConnectImp
 - CPositionImpDeviceNet, 321
- ConnectSlave
 - CGilsonDeviceNet, 97
- ControlState
 - HeadStageIDTypeState, 579
- CornerFrequency
 - CFilterPropertyNet, 71
- CornerFrequencymHz
 - CFilterPropertyNet, 71
- Count
 - CMcsUsbListNet, 223
- CPatchServerDeviceNet, 310
 - CPatchServerDeviceNet, 310
 - Sensor, 310
- CPathIdentDeviceNet, 311
 - ~CPathIdentDeviceNet, 311
 - CPathIdentDeviceNet, 311
 - Measure, 311
 - Set_Values, 311
- CPedoterDeviceNet, 312
 - !CPedoterDeviceNet, 312
 - ~CPedoterDeviceNet, 312
 - CPedoterDeviceNet, 312
 - GetCommand, 312
 - SetCommand, 313
- CPeristalticPumpDeviceNet, 313
 - ~CPeristalticPumpDeviceNet, 314
 - CPeristalticPumpDeviceNet, 314
 - McsBus_MotorControl, 314
- CPgaDeviceNet, 314
 - ~CPgaDeviceNet, 315
 - ApplyGains, 315
 - CPgaDeviceNet, 315
 - DefineAmplification, 315
 - DefineFrequencyRange, 315
 - DefineNumAmplifications, 315
 - DefineNumFrequencyRanges, 315
 - GetAmplification, 315
 - GetFrequencyRange, 316
 - GetGain, 316
 - GetNumAmplifications, 316
 - GetNumFrequencyRanges, 316
 - SetGain, 316
- CPositionIIDeviceNet, 316
 - !CPositionIIDeviceNet, 317
 - ~CPositionIIDeviceNet, 317
 - CPositionIIDeviceNet, 317
 - GetCoilCommunication, 317
 - GetImplantCurrentSetpoint, 318
 - GetImplantState, 318
 - GetOnOff, 318
 - RFFunction, 319
 - SetImplantCurrentSetpoint, 319
 - SwitchOnOff, 319
- CPositionImpDeviceNet, 319
 - !CPositionImpDeviceNet, 320
 - ~CPositionImpDeviceNet, 320
 - ConnectedImp, 321
 - ConnectImp, 321
 - CPositionImpDeviceNet, 320
 - GetDeviceList, 321
 - GetImpld, 321
 - GetRFFrequency, 321
 - SetDeviceList, 322
 - SetImpld, 322
 - SetRFFrequency, 322
- CPPCDeviceNet, 322
 - CPPCDeviceNet, 323
 - McsBus, 323
 - McsBus_MotorControl, 323
 - McsBus_Sensor, 323
 - PPCFunction, 323
- CPPCFunctionNet, 324
 - !CPPCFunctionNet, 325
 - ~CPPCFunctionNet, 325
 - CPPCFunctionNet, 325
 - FirePressurePulse, 325
 - GetAnalogVoltage, 326
 - GetAnalogVoltageRange, 326
 - GetDigitalIn, 326
 - GetPressureRange, 327
 - GetPumpModeType, 327
 - GetPumpSpeedUnit, 327
 - GetSupplyVoltage, 327
 - GetValveActive, 328
 - IsBusy, 328
 - LoadPressure, 328
 - MeasureReservoir, 329
 - SetAnalogVoltageRange, 329
 - SetPressureOffset, 329
 - SetPressureRange, 329
 - SetPumpModeType, 329
 - SetPumpSpeedUnit, 330
 - SetValveActive, 330
- CPPS_DeviceNet, 330

- CPPS_DeviceNet, 331
- McsBus, 331
- McsBus_MotorControl, 331
- McsBus_Sensor, 331
- PPS_Function, 331
- CPPS_FunctionNet, 331
 - CPPS_FunctionNet, 332
 - GetAnalogVoltage, 333
 - GetAnalogVoltages, 333
 - GetBubbleState, 333
 - GetDigitalIn, 333
 - GetPumpCouple, 333
 - GetPumpEnableSpeedRatio, 333
 - GetPumpFastOnOff, 333
 - GetPumpFastSpeed, 333
 - GetPumpFunctionSpeeds, 333
 - GetPumpManualOnOff, 333
 - GetPumpMaxSpeed, 334
 - GetPumpModeType, 334
 - GetPumpSpeedRatio, 334
 - GetPumpSpeedUnit, 334
 - GetSupplyVoltage, 334
 - GetUseBubble, 334
 - SetAnalogVoltages, 334
 - SetPumpCouple, 334
 - SetPumpEnableSpeedRatio, 334
 - SetPumpFastOnOff, 334
 - SetPumpFastSpeed, 335
 - SetPumpFunctionSpeeds, 335
 - SetPumpManualOnOff, 335
 - SetPumpMaxSpeed, 335
 - SetPumpModeType, 335
 - SetPumpSpeedRatio, 335
 - SetPumpSpeedUnit, 335
 - SetUseBubble, 335
- CPPSDeviceNet, 336
 - CPPSDeviceNet, 336
- CProgramPressureCurveNet, 336
 - !CProgramPressureCurveNet, 337
 - ~CProgramPressureCurveNet, 337
 - CProgramPressureCurveNet, 337
 - GetRepeats, 337
 - Program, 337
 - SetRepeats, 338
- CPulseGeneratorFunctionNet, 338
 - !CPulseGeneratorFunctionNet, 339
 - ~CPulseGeneratorFunctionNet, 339
 - CPulseGeneratorFunctionNet, 339
 - GetModeSelect, 339
 - GetPeriod, 340
 - GetPulseLength, 340
 - SetModeSelect, 340
 - SetPeriod, 340
 - SetPulseLength, 341
- CRadioControlledDevicesNet, 341
 - ConnectDevice, 342
 - CRadioControlledDevicesNet, 342
 - DisconnectDevice, 342
 - GetDeviceNames, 342
 - GetFrequency, 342
 - HasRadioControl, 342
 - SetFrequency, 342
 - StillConnected, 342
- CreateSideband
 - CStimulusFunctionNet, 462
- CreateWirelessHeadstageSerialNumberString
 - CWirelessBaseFunctionNet, 563
- CRegionOfInterestRect
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, 343
- CRetinaLedDeviceNet, 344
 - ~CRetinaLedDeviceNet, 344
 - AddLoopEntry, 344
 - AddTableEntry, 344
 - ClearTable, 345
 - CRetinaLedDeviceNet, 344
 - GetTablepointer, 345
 - SetLED, 345
 - SetLumi, 345
 - SetPersistency, 345
 - SetRepeat, 345
 - SetTablepointer, 345
 - SetTrigger, 345
- CRFFFunctionNet, 346
 - !CRFFFunctionNet, 347
 - ~CRFFFunctionNet, 347
 - Connect, 347
 - CRFFFunctionNet, 347
 - GetAvailableDeviceList, 347
 - GetAvailableDeviceListEx, 347
 - GetAvailableStateList, 348
 - GetAvailableStateListEx, 348
 - GetBaseFrequency, 348
 - GetConnectedDevice, 349
 - GetState, 349
 - GetTestMode, 349
 - GetWorkingFrequency, 349
 - SetBaseFrequency, 349
 - SetTestMode, 350
 - SetWorkingFrequency, 350
- CRobo_FYIProgram_FunctionNet, 350
 - CRobo_FYIProgram_FunctionNet, 351
 - GetLength, 351
 - GetState, 351
 - GetValve1, 351
 - GetValve2, 351
 - SetLength, 351
 - SetValve1, 351
 - SetValve2, 352
 - Start, 352
- CRobo_FYITemp_FunctionNet, 352
 - CRobo_FYITemp_FunctionNet, 352
 - GetICoeff, 353
 - GetMaxPower, 353
 - GetPCoeff, 353
 - GetRegulatorOnOff, 353

- GetSollTemp, [353](#)
- SetICoeff, [353](#)
- SetMaxPower, [353](#)
- SetPCoeff, [353](#)
- SetRegulatorOnOff, [353](#)
- SetSollTemp, [354](#)
- CRoboDacqNet, [354](#)
 - CancelTableLoop, [356](#)
 - CancelTableLoopAndStopTable, [356](#)
 - ClampAmpRestart, [357](#)
 - CRoboDacqNet, [356](#)
 - DoRamp, [357](#)
 - Emu_GetCellCapacity, [357](#)
 - Emu_GetCellPotential, [357](#)
 - Emu_GetCellResists, [357](#)
 - Emu_GetElectrodeResists, [357](#)
 - Emu_GetNoise, [357](#)
 - Emu_SetCellCapacity, [357](#)
 - Emu_SetCellPotential, [357](#)
 - Emu_SetCellResists, [357](#)
 - Emu_SetElectrodeResists, [358](#)
 - Emu_SetNoise, [358](#)
 - GetAllDigout, [358](#)
 - GetCapacityC, [358](#)
 - GetCapacityV, [358](#)
 - GetCapacityX, [358](#)
 - GetClampAmpSerialNumber, [358](#)
 - GetCommand, [358](#)
 - GetConfigurationBit, [358](#)
 - GetConfigurationBitAxc, [358](#)
 - GetConfigurationBitBlu_Led, [359](#)
 - GetConfigurationBitBlu_LedToggleFast, [359](#)
 - GetConfigurationBitBlu_LedToggleSlow, [359](#)
 - GetConfigurationBitCC_Gen, [359](#)
 - GetConfigurationBitCV_Gen, [359](#)
 - GetConfigurationBitRC_Gen, [359](#)
 - GetConfigurationBitRed_Led, [359](#)
 - GetConfigurationBitRed_LedSaturation, [359](#)
 - GetConfigurationBitRed_LedToggleFast, [359](#)
 - GetConfigurationBitRed_LedToggleSlow, [359](#)
 - GetConfigurationBitRelais, [359](#)
 - GetConfigurationBitRV_Gen, [360](#)
 - GetConfigurationBits, [360](#)
 - GetConfigurationBitStream, [360](#)
 - GetConfigurationBitSupply, [360](#)
 - GetCrossTalkOffset, [360](#)
 - GetCrossTalkOptimum, [360](#)
 - GetDigout, [360](#)
 - GetDisplayText, [360](#)
 - GetDownsampleFactor, [360](#)
 - GetFilter, [360](#)
 - GetFilterCoeffs, [360](#)
 - GetIC, [361](#)
 - GetIClamp, [361](#)
 - GetICOffset, [361](#)
 - GetIGain, [361](#)
 - GetNIC_MS, [361](#)
 - GetNUC_MS, [361](#)
- GetNUV_MS, [361](#)
- GetPGain, [361](#)
- GetRecordingNumber, [361](#)
- GetResistanceC, [361](#)
- GetResistanceV, [361](#)
- GetScreen, [362](#)
- GetSimulation, [362](#)
- GetUC, [362](#)
- GetUClamp, [362](#)
- GetUCOffset, [362](#)
- GetUpdateDisplay, [362](#)
- GetUV, [362](#)
- GetUVOffset, [362](#)
- GetXGain, [362](#)
- RunTable, [362](#)
- SetAllDigout, [363](#)
- SetCommand, [363](#)
- SetConfigurationBit, [363](#)
- SetConfigurationBitAxc, [363](#)
- SetConfigurationBitBlu_Led, [363](#)
- SetConfigurationBitBlu_LedToggleFast, [363](#)
- SetConfigurationBitBlu_LedToggleSlow, [363](#)
- SetConfigurationBitCC_Gen, [363](#)
- SetConfigurationBitCV_Gen, [363](#)
- SetConfigurationBitRC_Gen, [364](#)
- SetConfigurationBitRed_Led, [364](#)
- SetConfigurationBitRed_LedSaturation, [364](#)
- SetConfigurationBitRed_LedToggleFast, [364](#)
- SetConfigurationBitRed_LedToggleSlow, [364](#)
- SetConfigurationBitRelais, [364](#)
- SetConfigurationBitRV_Gen, [364](#)
- SetConfigurationBitStream, [364](#)
- SetConfigurationBitSupply, [364](#)
- SetCrossTalkOffset, [364](#)
- SetCrossTalkOptimum, [365](#)
- SetDigout, [365](#)
- SetDisplayText, [365](#)
- SetDownsampleFactor, [365](#)
- SetFilter, [365](#)
- SetFilterCoeffs, [365](#)
- SetIClamp, [365](#)
- SetICOffset, [365](#)
- SetIGain, [365](#)
- SetPGain, [366](#)
- SetRecordingNumber, [366](#)
- SetScreen, [366](#)
- SetSimulation, [366](#)
- SetUClamp, [366](#)
- SetUCOffset, [366](#)
- SetUVOffset, [366](#)
- SetXGain, [366](#)
- StopTable, [366](#)
- Table_Wait, [367](#)
- TableDefBegin, [367](#)
- TableDefEnd, [367](#)
- UpdateDisplay, [367](#)
- CRoboDeviceNet, [367](#)
 - ~CRoboDeviceNet, [370](#)

- Axes_I, [377](#)
- Axes_X, [377](#)
- Axes_Y, [377](#)
- Axes_Z, [377](#)
- Axis_I, [377](#)
- Axis_X, [377](#)
- Axis_Y, [377](#)
- Axis_Z, [378](#)
- CancelPoolLoop, [371](#)
- CancelPoolLoopAndStopMovement, [371](#)
- CRoboDeviceNet, [370](#)
- FindReference, [371](#)
- GetAirpressure, [371](#)
- GetAirpressureLimit, [371](#)
- GetAirValve, [371](#)
- GetCurrentAirvalve, [371](#)
- GetCurrentAirvalveLimit, [372](#)
- GetCurrentPosition, [372](#)
- GetErrorAirpressure, [372](#)
- GetErrorCurrentAirvalve, [372](#)
- GetErrorVoltage12V, [372](#)
- GetErrorVoltage5V, [372](#)
- GetErrorVoltageAirvalve, [372](#)
- GetErrorVoltageRs485A, [372](#)
- GetErrorVoltageRs485B, [373](#)
- GetErrorVoltageValves, [373](#)
- GetInMovement, [373](#)
- GetMinPressure, [373](#)
- GetMovementError, [373](#)
- GetVoltage12V, [373](#)
- GetVoltage12VLimit, [373](#)
- GetVoltage5V, [373](#)
- GetVoltage5VLimit, [373](#)
- GetVoltageAirvalve, [373](#)
- GetVoltageAirvalveLimit, [374](#)
- GetVoltageRs485A, [374](#)
- GetVoltageRs485ALimit, [374](#)
- GetVoltageRs485B, [374](#)
- GetVoltageRs485BLimit, [374](#)
- GetVoltageValves, [374](#)
- GetVoltageValvesLimit, [374](#)
- McsBus, [381](#)
- McsBus_MotorControl, [381](#)
- McsBus_XY, [378](#)
- McsBus_ZI, [378](#)
- MoveAbs, [374](#)
- RoboError_AnotherMaster, [378](#)
- RoboError_Base, [378](#)
- RoboError_CannotEscapeEndSwitch, [378](#)
- RoboError_CommandAlreadyInProgress, [378](#)
- RoboError_CommandNotPossible, [378](#)
- RoboError_CommunicationTimeout, [378](#)
- RoboError_DacqNotReady, [379](#)
- RoboError_DLLMovementTimeout, [379](#)
- RoboError_FindReferenceMethod, [379](#)
- RoboError_GilsonCommandPending, [379](#)
- RoboError_GilsonTimeout, [379](#)
- RoboError_GilsonWronId, [379](#)
- RoboError_McsBus_UnknownCommand, [379](#)
- RoboError_NoEndSwitch, [379](#)
- RoboError_NoMoreData, [379](#)
- RoboError_NoReference, [379](#)
- RoboError_NoSpeedOrAcceleration, [380](#)
- RoboError_OverPressure, [380](#)
- RoboError_ParameterNotAllowed, [380](#)
- RoboError_PeristalticTimeout, [380](#)
- RoboError_Phase0OutOfRange, [380](#)
- RoboError_PollLoopCanceled, [380](#)
- RoboError_PollLoopCanceledAndStopMovement, [380](#)
- RoboError_Pressure, [380](#)
- RoboError_RangeExceeded, [380](#)
- RoboError_StateChangeNotPossible, [380](#)
- RoboError_Timeout, [381](#)
- RoboError_UnknownCommand, [381](#)
- RoboMainLowLevelCommand, [381](#)
- RoboStatusEvent, [381](#)
- SetAirpressureLimit, [375](#)
- SetAirValve, [375](#)
- SetCurrentAirvalveLimit, [375](#)
- SetCurrentAndAir, [375](#)
- SetInMovement, [375](#)
- SetMinPressure, [375](#)
- SetVoltage12VLimit, [376](#)
- SetVoltage5VLimit, [376](#)
- SetVoltageAirvalveLimit, [376](#)
- SetVoltageRs485ALimit, [376](#)
- SetVoltageRs485BLimit, [376](#)
- SetVoltageValvesLimit, [376](#)
- StopMovement, [376](#)
- CRoboDeviceNet::RoboMainLowLevelCommands, [583](#)
- FindReferencePhase0, [584](#)
- GetAxisConfig, [584](#)
- GetHWConfig, [584](#)
- GetHWRevision, [584](#)
- GetMaxNoPressure, [584](#)
- GetMaxNoPressureWaitTime, [584](#)
- GetMaxPressureWaitTime, [584](#)
- GetMinNoPressureWaitTime, [584](#)
- GetMinPressure, [584](#)
- GetMinPressureWaitTime, [585](#)
- GetParameter, [585](#)
- GetPhases, [585](#)
- GetSearchReferenceFastAccel, [585](#)
- GetSearchReferenceFastSpeed, [585](#)
- GetSearchReferenceFineAccel, [585](#)
- GetSearchReferenceFineSpeed, [585](#)
- GetSearchReferenceMethod, [586](#)
- GetSearchReferenceMoveOut, [586](#)
- GetSearchReferenceOffsetPos, [586](#)
- GetUserParameter, [586](#)
- HasRef, [587](#)
- SetAxisConfig, [587](#)
- SetHWConfig, [587](#)
- SetHWRevision, [587](#)
- SetMaxNoPressure, [587](#)

- SetMaxNoPressureWaitTime, [587](#)
- SetMaxPressureWaitTime, [587](#)
- SetMinNoPressureWaitTime, [587](#)
- SetMinPressure, [587](#)
- SetMinPressureWaitTime, [587](#)
- SetParameter, [588](#)
- SetSearchReferenceFastAccel, [588](#)
- SetSearchReferenceFastSpeed, [588](#)
- SetSearchReferenceFineAccel, [588](#)
- SetSearchReferenceFineSpeed, [588](#)
- SetSearchReferenceMethod, [588](#)
- SetSearchReferenceMoveOut, [589](#)
- SetSearchReferenceOffsetPos, [589](#)
- SetUserParameter, [589](#)
- CRoboFluidDeviceNet, [381](#)
 - ~CRoboFluidDeviceNet, [382](#)
 - CloseAllValves, [382](#)
 - CRoboFluidDeviceNet, [382](#)
 - GetPumpSpeed, [383](#)
 - GetSingleValve, [383](#)
 - GetValve, [383](#)
 - IsPumpMotorOn, [383](#)
 - m_pMcsBus_MotorControlNet, [384](#)
 - m_pRoboFluidDevice, [384](#)
 - McsBus_MotorControl, [384](#)
 - PumpOff, [383](#)
 - PumpOn, [383](#)
 - SetPumpSpeed, [383](#)
 - SetSingleValve, [384](#)
 - SetValve, [384](#)
- CRoboInjectDeviceNet, [385](#)
 - CRoboInjectDeviceNet, [385](#)
- CRoboocyte2DeviceNet, [385](#)
 - CRoboocyte2DeviceNet, [386](#)
 - GetAxisLED, [386](#)
 - GetGilsonDevice, [386](#)
 - GetMcsBus_Extension, [386](#)
 - GetRoboDacq, [386](#)
 - GetRoboFluidDevice, [386](#)
 - SetAxisLED, [387](#)
- CRoboStatorDeviceNet, [387](#)
 - CRoboStatorDeviceNet, [388](#)
 - FindReferenceI, [388](#)
 - FindReferenceXY, [388](#), [389](#)
 - FindReferenceZ, [389](#)
 - GetCurrentPositionI, [389](#)
 - GetCurrentPositionXY, [389](#)
 - GetCurrentPositionZ, [389](#)
 - HasRefI, [389](#)
 - HasRefXY, [389](#)
 - HasRefZ, [389](#)
 - MoveAbsI, [389](#), [390](#)
 - MoveAbsXY, [390](#)
 - MoveAbsZ, [390](#)
 - RoboMainStatorLowLevelCommand, [393](#)
 - SetAccelerationI, [390](#)
 - SetAccelerationNativeI, [390](#)
 - SetAccelerationNativeXY, [390](#)
 - SetAccelerationNativeZ, [390](#)
 - SetAccelerationXY, [391](#)
 - SetAccelerationZ, [391](#)
 - SetCurrentAndAirXY, [391](#)
 - SetSpeedI, [391](#)
 - SetSpeedNativeI, [391](#)
 - SetSpeedNativeXY, [391](#)
 - SetSpeedNativeZ, [391](#)
 - SetSpeedXY, [391](#)
 - SetSpeedZ, [392](#)
 - SetVelocityI, [392](#)
 - SetVelocityXY, [392](#)
 - SetVelocityZ, [392](#)
 - StopMovementI, [392](#)
 - StopMovementXY, [392](#)
 - StopMovementZ, [392](#)
- CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands, [590](#)
 - FindReferencePhase0XY, [590](#)
- CSafeISDeviceNet, [393](#)
 - ~CSafeISDeviceNet, [394](#)
 - CSafeISDeviceNet, [394](#)
 - DacqDevice, [395](#)
 - FluidControlDevice, [396](#)
 - RoboDevice, [396](#)
 - SetAdcChannels, [394](#)
 - SetAdcSamplePos, [394](#)
 - SetDacMode, [394](#)
 - SetDacPeriode, [395](#)
 - SetDacPulseform, [395](#)
 - SetSwitches, [395](#)
- CSCUDacqGroupChannelSelectionNet, [396](#)
 - CSCUDacqGroupChannelSelectionNet, [396](#)
- CSCUFunctionNet, [397](#)
 - !CSCUFunctionNet, [399](#)
 - ~CSCUFunctionNet, [399](#)
 - AutomaticAnalogOut, [399](#)
 - CSCUFunctionNet, [399](#)
 - EnableAnalogOut, [400](#)
 - GetAnalogOutADCRange, [400](#)
 - GetAnalogOutChannels, [400](#)
 - GetAnalogOutDACRange, [400](#)
 - GetAvailableHeadstages, [400](#)
 - GetAvailableHeadstagesEvent, [412](#)
 - GetFilterProperties, [401](#)
 - GetFilterProperty, [401](#)
 - GetHeadstageAdcBits, [401](#)
 - GetHeadstageAdcRangeInMicroVolt, [402](#)
 - GetHeadstageDacBits, [402](#)
 - GetHeadstageDacCurrentRangeInMicroAmpere, [402](#)
 - GetHeadstageDacCurrentResolutionInNanoAmpere, [403](#)
 - GetHeadstageDacVoltageRangeInMilliVolt, [403](#)
 - GetHeadstageDacVoltageResolutionInMicroVolt, [403](#)
 - GetHeadstageGainInPer mille, [404](#)
 - GetHeadstageID, [404](#)

- GetHeadstageNumberOfAnalogChannels, 404
- GetHeadstageNumberOfStimulationChannels, 405
- GetHeadstageSamplerate, 405
- GetHeadstageSerialNumber, 405
- GetMaxNumberOfHeadstages, 406
- GetMaxStimulusChannelsPerHeadstage, 406
- GetReferenceElectrodeMode, 406
- GetReferenceElectrodeSwitchState, 406
- HasAnalogOut, 407
- HasGalvanicIsolation, 407
- HasHSPowerSwitch, 407
- IsAnalogOutEnabled, 407
- IsAutomaticAnalogOut, 408
- IsHeadstageAvailable, 408
- IsHeadstageAvailableEvent, 412
- IsHSPowered, 408
- IsInDacqLegacyMode, 408
- OnGetAvailableHeadstages, 409
- OnIsHeadstageAvailable, 409
- PowerHS, 409
- SetAnalogOutADCRange, 409
- SetAnalogOutChannels, 409
- SetAnalogOutDACRange, 411
- SetDacqLegacyMode, 411
- SetReferenceElectrodeMode, 411
- SetReferenceElectrodeSwitchState, 411
- CSerialPortNet, 412
 - CSerialPortNet, 412
 - GetBytesAvailable, 413
 - Receive, 413
 - ReceiveString, 413
 - Send, 413
- CStg200xBasicNet, 414
 - ~CStg200xBasicNet, 418
 - GetAnalogRanges, 418
 - GetAnalogResolution, 418
 - GetAutocalibrationDisabled, 419
 - GetAvailableMemory, 419
 - GetBlankingEnable, 419, 420
 - GetCurrentRangeInNanoAmp, 420
 - GetCurrentResolutionInNanoAmp, 420
 - GetDacAmplificationFactor, 421
 - GetDACResolution, 421
 - GetDiginValue, 421
 - GetDigoutMode, 421
 - GetDigoutValue, 422
 - GetElectrodeDacMux, 422
 - GetElectrodeEnable, 423
 - GetElectrodeMode, 424
 - GetEnableAmplifierProtectionSwitch, 424, 425
 - GetExternalElectrodeEnable, 425
 - GetFAAmplification, 426
 - GetHeadstage, 426
 - GetListModelIndexRange, 426
 - GetListmodeTriggerSource, 426
 - GetNumberOfAnalogChannels, 426
 - GetNumberOfHWDACPaths, 426
 - GetNumberOfStimulationElectrodes, 426
 - GetNumberOfStimulationSourcesPerElectrode, 427
 - GetNumberOfSyncoutChannels, 427
 - GetNumberOfTriggerInputs, 427
 - GetOutputRate, 427
 - GetStgProgramInfo, 427, 428
 - GetStgVersionInfo, 428
 - GetSyncoutMap, 428
 - GetTotalMemory, 429
 - GetTriggerSource, 429
 - GetVoltageRangeInMicroVolt, 429
 - GetVoltageResolutionInMicroVolt, 429
 - ListModeSendStart, 430
 - ListModeSendStop, 430
 - SendStart, 430
 - SendStop, 430
 - SetAutocalibrationDisabled, 430
 - SetBlankingEnable, 431, 432
 - SetCurrentMode, 432
 - SetDacAmplificationFactor, 432
 - SetDigoutMode, 433
 - SetDigoutValue, 433
 - SetElectrodeDacMux, 433–435
 - SetElectrodeEnable, 435–437
 - SetElectrodeMode, 438, 439
 - SetEnableAmplifierProtectionSwitch, 439, 440
 - SetExternalElectrodeEnable, 441
 - SetFAAmplification, 442
 - SetHeadstage, 442
 - SetListModelIndexRange, 442
 - SetListmodeTriggerSource, 442
 - SetMeasurementMode, 442
 - SetOutputRate, 443
 - SetStgProgramInfo, 443
 - SetSyncoutMap, 443
 - SetTriggerSource, 443, 444
 - SetVoltageMode, 444
- CStg200xDownloadBasicNet, 444
 - ClearChannelData, 446
 - ClearSyncData, 446
 - DisableAutoReset, 446
 - EnableAutoReset, 447
 - ForceStatusEvent, 447
 - GetMemoryUsageDAC, 447
 - GetMemoryUsageSyncout, 447
 - GetSweepCount, 448
 - GetTrigger, 448
 - ResetStatus, 448
 - SendChannelData, 450
 - SendSyncData, 450
 - SetupRetriggerMode, 451
 - SetupTrigger, 451
 - SetupTriggerSingle, 452
 - Stimulus, 453
- CStg200xDownloadNet, 453
 - ~CStg200xDownloadNet, 454
 - ClearChannel_PrepAndSendData, 455
 - CStg200xDownloadNet, 454

- DisableMultiFileMode, [455](#)
- EnableMultiFileMode, [455](#)
- GetModuleCurrent, [456](#)
- GetModuleTemp, [456](#)
- MwPollStatusEvent, [459](#)
- PrepareAndAppendData, [456](#)
- PrepareAndSendData, [457](#)
- QueryTriggerstatus, [458](#)
- SendSegmentDefine, [458](#)
- SendSegmentSelect, [458](#)
- SendSegmentStart, [459](#)
- SetOutputMap, [459](#)
- Stg200xPollStatusEvent, [459](#)
- CStimulusFunctionNet, [460](#)
 - ClearChannel_PrepareAndSendData, [462](#)
 - ClearChannelData, [462](#)
 - ClearMultiplexedData, [462](#)
 - ClearSyncData, [462](#)
 - CreateSideband, [462](#)
 - CStimulusFunctionNet, [461](#)
 - ForceStatusEvent, [463](#)
 - GetAvailableMemory, [463](#)
 - GetCurrentRangeInNanoAmp, [464](#)
 - GetCurrentResolutionInNanoAmp, [464](#)
 - GetDACResolution, [464](#)
 - GetMultiplexedDataChannelsInBlock, [464](#)
 - GetNumberOfAnalogChannels, [465](#)
 - GetTotalMemory, [465](#)
 - GetVoltageRangeInMicroVolt, [465](#)
 - GetVoltageResolutionInMicroVolt, [465](#)
 - PollStatusEvent, [471](#)
 - PrepareAndAppendData, [466](#)
 - PrepareAndSendData, [467](#)
 - PrepareData, [467](#)
 - SendMultiplexedData, [468](#)
 - SendPreparedData, [468](#)
 - SendStart, [468](#)
 - SendStop, [468](#)
 - SetupTrigger, [469](#)
 - SetupTriggerSingle, [469](#)
 - StartPoll, [471](#)
 - StopPoll, [471](#)
- CStimulusFunctionNet::SidebandData, [590](#)
 - !SidebandData, [591](#)
 - ~SidebandData, [591](#)
 - Duration, [591](#)
 - Sideband, [591](#)
 - SidebandData, [591](#)
- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [592](#)
 - !StimulusDeviceDataAndUnrolledData, [593](#)
 - ~StimulusDeviceDataAndUnrolledData, [593](#)
 - DeviceData, [593](#)
 - DeviceDataLength, [593](#)
 - StimulusDeviceDataAndUnrolledData, [592](#)
 - UnrolledAmplitude, [593](#)
 - UnrolledDuration, [593](#)
 - UnrolledSync, [593](#)
- CSw2to64DeviceNet, [471](#)
 - ~CSw2to64DeviceNet, [472](#)
 - CSw2to64DeviceNet, [472](#)
 - GetChannel, [472](#)
 - GetChannels, [473](#)
 - GetNumber, [473](#)
 - SetChannel, [473](#)
 - SetChannels, [473](#)
- CTcxDeviceNet, [474](#)
 - ~CTcxDeviceNet, [476](#)
 - CalibrateThermocouple, [476](#)
 - CTcxDeviceNet, [476](#)
 - FactoryReset, [476](#)
 - GetBoardTemp, [477](#)
 - GetCalibration, [477](#)
 - GetCalibrationDecp, [477](#)
 - GetCalibrationMax, [477](#)
 - GetCalibrationMin, [477](#)
 - GetCurrent, [477](#)
 - GetD, [477](#)
 - GetDDecp, [477](#)
 - GetDevice, [477](#)
 - GetDeviceType, [478](#)
 - GetDevname, [478](#)
 - GetDMax, [478](#)
 - GetDMin, [478](#)
 - GetDuty, [478](#)
 - GetEnableHeaterLimit, [478](#)
 - GetEnableThermocouple, [478](#)
 - GetHasThermocouple, [478](#)
 - GetHeaterLimit, [479](#)
 - GetHeaterTemp, [479](#)
 - GetI, [479](#)
 - GetIDecp, [479](#)
 - GetIMax, [479](#)
 - GetIMin, [479](#)
 - GetIOut, [479](#)
 - GetMaxHeaterPowerMultiwell, [480](#)
 - GetMaxP, [480](#)
 - GetMaxpDecp, [480](#)
 - GetMaxpMax, [480](#)
 - GetMaxpMin, [480](#)
 - GetNumControlChannels, [480](#)
 - GetNumDevices, [480](#)
 - GetNumMeasureChannels, [480](#)
 - GetOnOff, [481](#)
 - GetP, [481](#)
 - GetPDecp, [481](#)
 - GetPMax, [481](#)
 - GetPMin, [481](#)
 - GetPOut, [481](#)
 - GetPwrOut, [481](#)
 - GetPwrSet, [482](#)
 - GetRes1, [482](#)
 - GetRes2, [482](#)
 - GetResS, [482](#)
 - GetResX, [482](#)
 - GetROut, [482](#)

- GetSensorType, [482](#)
- GetSetpoint, [482](#)
- GetSetpointDecp, [483](#)
- GetSetpointMax, [483](#)
- GetSetpointMin, [483](#)
- GetThermocoupleCalibration, [483](#)
- GetThermocoupleNanovoltPerKelvin, [483](#)
- GetThermocoupleReferenceTemp, [483](#)
- GetThermocoupleTemp, [483](#)
- GetThermocoupleTempAbs, [484](#)
- GetUnit, [484](#)
- GetUOut, [484](#)
- GetValue, [484](#)
- GetValueHires, [484](#)
- GetVolti, [484](#)
- SetCalibration, [484](#)
- SetD, [485](#)
- SetDevice, [485](#)
- SetDeviceType, [485](#)
- SetDevname, [485](#)
- SetEnableHeaterLimit, [485](#)
- SetEnableThermocouple, [485](#)
- SetHeaterLimit, [485](#)
- SetI, [485](#)
- SetMaxHeaterPowerMultiwell, [486](#)
- SetMaxP, [486](#)
- SetOnOff, [486](#)
- SetP, [486](#)
- SetSensorType, [486](#)
- SetSetpoint, [486](#)
- SetThermocoupleNanovoltPerKelvin, [487](#)
- CTEERFunctionNet, [487](#)
 - !CTEERFunctionNet, [489](#)
 - ~CTEERFunctionNet, [489](#)
- CancelInternalCalibration, [489](#)
- CTEERFunctionNet, [489](#)
- GetAdapterCode, [489](#)
- GetAdcOffsetU1, [490](#)
- GetAdcOffsetU2, [490](#)
- GetAmplitude_nA, [490](#)
- GetBytesPerSample, [490](#)
- GetClampMode, [490](#)
- GetControllerParams, [490](#)
- GetCurrentEnable, [491](#)
- GetDacZero, [491](#)
- GetLiquidResistance, [491](#)
- GetMaxChunkSize_Byte, [491](#)
- GetNumberOfAvailableSamples, [491](#)
- GetPeriod_us, [492](#)
- GetRotaryPositionCode, [492](#)
- GetSampleBufferChunk, [492](#)
- GetSampleVoltageBuffer_uV, [492](#)
- GetScaleFactorU1, [493](#)
- GetScaleFactorU2, [493](#)
- GetUptimeSeconds, [493](#)
- GetWaveform, [493](#)
- IsInternalCalibrationFinished, [493](#)
- IsSamplingFinished, [494](#)
- SetAmplitude_nA, [494](#)
- SetBufferIndex, [494](#)
- SetClampMode, [494](#)
- SetControllerParams, [495](#)
- SetCurrentEnable, [495](#)
- SetExternalLED, [495](#)
- SetLiquidResistance, [495](#)
- SetPeriod_us, [496](#)
- SetWaveform, [496](#)
- StartInternalCalibration, [496](#)
- StartSampling, [496](#)
- StopSampling, [496](#)
- CTEERMachineDeviceNet, [497](#)
 - ~CTEERMachineDeviceNet, [497](#)
 - CTEERMachineDeviceNet, [497](#)
 - TEERFunctionNet, [497](#)
- CurrentRangeInNanoAmp
 - W2100_StimulusParametersNet, [594](#)
- CurrentResolutionInNanoAmp
 - W2100_StimulusParametersNet, [594](#)
- CUsbExceptionNet, [498](#)
 - CUsbExceptionNet, [498](#)
 - Status, [499](#)
- CutoffFrequency
 - CCreateFilterNet, [47](#)
- CVoltageRangeInfoNet
 - CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, [499](#)
- CW2100_FunctionNet, [499](#)
 - ClearStimulusParametersCache, [501](#)
 - ClearUserDefinedNameCache, [501](#)
 - CW2100_FunctionNet, [501](#)
 - DeselectAllHeadstages, [502](#)
 - DeselectHeadstage, [502](#)
 - GetAccelGyroCurrentRate, [502](#)
 - GetAccelGyroDesiredRate, [502](#)
 - GetAccelGyroEnabled, [502](#)
 - GetAccelRange, [502](#)
 - GetAnalogOutChannel, [502](#)
 - GetAnalogOutFilter, [502](#)
 - GetAudioChannels, [502](#)
 - GetAvailableHeadstages, [502](#)
 - GetBatteryState, [503](#)
 - GetDacRange, [503](#)
 - GetFilterProperties, [503](#)
 - GetFilterProperty, [503](#)
 - GetFPGA FirmwareType, [503](#)
 - GetGyroRange, [503](#)
 - GetHeadstageOnOff, [503](#)
 - GetHeadstageSamplingActive, [503](#)
 - GetMultiHeadstageMode, [503](#)
 - GetPicFirmwareType, [503](#)
 - GetSelectedChannels, [504](#)
 - GetSelectedHeadstageState, [504](#)
 - GetStimulusParametersCache, [504](#)
 - GetStimulusParametersFromSelectedHS, [504](#)
 - GetStiumlusParameters, [504](#)
 - GetUserDefinedName, [504](#)

- GetUserDefinedNameCache, 504
- GetUserDefinedNameFromSelectedHS, 505
- PulseGenerator, 506
- SelectHeadstage, 505
- SetAccelGyroDesiredRate, 505
- SetAccelGyroEnabled, 505
- SetAccelRange, 505
- SetAnalogOutChannel, 505
- SetAnalogOutFilter, 505
- SetAudioChannels, 505
- SetDacRange, 506
- SetGyroRange, 506
- SetHeadstageOnOff, 506
- SetHeadstageSamplingActive, 506
- SetHeadstageToSleep, 506
- SetMultiHeadstageMode, 506
- SetSelectedChannels, 506
- Stimulator, 507
- CW2100_FunctionNet::AudioChannelsNet, 27
 - amplification, 27
 - channel, 27
 - dacqgroup, 27
- CW2100_StimulatorFunctionNet, 507
 - BOOST_BIT, 513
 - ClearChannelData, 508
 - CW2100_StimulatorFunctionNet, 508
 - GetBoostAlwaysOnMode, 509
 - GetBoostPreTime, 509
 - GetCurrentRangeInNanoAmp, 509
 - GetCurrentResolutionInNanoAmp, 509
 - GetDACResolution, 509
 - GetDigitalStimulatorTrigger, 509
 - GetDigitalStimulatorTriggerSlope, 510
 - GetNumberOfAnalogChannels, 510
 - GetNumberOfSyncoutChannels, 510
 - GetNumberOfTriggerInputs, 510
 - GetStimulationPatternMemory, 510
 - GetTimeResolutionInNanoSeconds, 510
 - GetTimeSlot, 510
 - GetVoltageRangeInMicroVolt, 510
 - GetVoltageResolutionInMicroVolt, 511
 - GND_SWITCH_BIT, 513
 - PollStatusEvent, 513
 - PrepareData, 511
 - PrepareDataSync, 511
 - SelectTimeSlot, 511
 - SendPreparedData, 511
 - SendStart, 512
 - SendStop, 512
 - SetBoostAlwaysOnMode, 512
 - SetDigitalStimulatorTrigger, 512
 - SetDigitalStimulatorTriggerSlope, 512
 - StartPoll, 512
 - StopPoll, 513
 - SYNC_BIT0, 513
 - SYNC_BIT1, 513
- CW2100DacqGroupChannelSelectionNet, 513
 - CW2100DacqGroupChannelSelectionNet, 514
- CWarnerUssingDeviceNet, 514
 - !CWarnerUssingDeviceNet, 515
 - ~CWarnerUssingDeviceNet, 515
 - CWarnerUssingDeviceNet, 515
 - WarnerUssingFunction, 515
- CWarnerUssingFunctionNet, 515
 - !CWarnerUssingFunctionNet, 518
 - ~CWarnerUssingFunctionNet, 518
 - CompensateElectrodeOffset, 518
 - CWarnerUssingFunctionNet, 517
 - GetAvailableChambers, 518
 - GetChannelsCountOfChamber, 518
 - GetClampMode, 519
 - GetComplianceVoltage, 519
 - GetDacNampsPerDigitHighCurrentRange, 519
 - GetDacPampsPerDigitLowCurrentRange, 520
 - GetDacZero, 520
 - GetHighCurrentRange, 520
 - GetIdleModeOffset, 521
 - GetLiquidResistance, 521
 - GetLowCurrentRange, 521
 - GetNumberOfAvailableChambers, 522
 - GetNumberOfHardwareSlotsForChambers, 522
 - GetU1Offset, 522
 - GetU1Reference, 522
 - GetU2Offset, 523
 - GetU2Reference, 523
 - GetUnitDescription, 523
 - GetUnitExponent, 524
 - GetUnitName, 524
 - GetUnitsPerDigit, 524
 - GetUptimeSeconds, 526
 - GetVoltageClampControllerParam_D, 526
 - GetVoltageClampControllerParam_I, 526
 - GetVoltageClampControllerParam_P, 527
 - IsChamberAvailable, 527
 - IsHighCurrentMode, 527
 - IsInternalCalibrationFinished, 528
 - IsPulseEnabled, 528
 - SetClampMode, 528
 - SetEnablePulse, 529
 - SetHighCurrentMode, 529
 - SetIdleModeOffset, 529
 - SetLiquidResistance, 529
 - SetLowCurrentMode, 530
 - SetPulse, 530
 - SetVoltageClampControllerParam_D, 530
 - SetVoltageClampControllerParam_I, 531
 - SetVoltageClampControllerParam_P, 531
 - WaitForAllChambers, 531
 - WaitForChamber, 531
- CWarnerValveControllerDeviceNet, 532
 - !CWarnerValveControllerDeviceNet, 536
 - ~CWarnerValveControllerDeviceNet, 536
 - ClearTableName, 536
 - ClearValveTable, 536
 - CWarnerValveControllerDeviceNet, 536
 - GetActiveRunningTableNumber, 536

- GetActiveRunningTableNumberEvent, 552
- GetAnalogThresholdHigh, 537
- GetAnalogThresholdHighEvent, 552
- GetAnalogThresholdLow, 537
- GetAnalogThresholdLowEvent, 552
- GetAnalogVoltage, 537
- GetAnalogVoltageEvent, 552
- GetCurrentEditTableNumber, 537
- GetCurrentNumberOfValves, 538
- GetCurrentNumberOfValvesEvent, 552
- GetDigitalOutPortValve, 538
- GetDigitalOutPortValveEvent, 553
- GetDigitalPortDirection, 538
- GetDigitalPortDirectionEvent, 553
- GetDisplayMode, 538
- GetDisplayModeEvent, 553
- GetTableName, 539
- GetTableNamebyIndex, 539
- GetTableNamebyIndexEvent, 553
- GetTotalNumberOfDigitalPorts, 539
- GetTotalNumberOfTables, 539
- GetTotalNumberOfValves, 539
- GetTotalTableSize, 540
- GetValveActive, 540
- GetValveActiveEvent, 553
- GetValveBoardRevision, 540
- GetValveBoardRevisionEvent, 553
- GetValveBoardRevisionString, 540
- GetValveDigitalInPort, 540
- GetValveDigitalInPortEvent, 553
- GetValveManualGroup, 541
- GetValveManualGroupEvent, 554
- GetValveManualState, 541
- GetValveManualStateEvent, 554
- GetValveMode, 541
- GetValveModeEvent, 554
- GetValvesActiveMap, 542
- GetValvesManualStateMap, 542
- GetValveTableEntry, 542
- IsDigitalOutPortInverted, 542
- IsDigitalOutPortInvertedEvent, 554
- IsValveDigitalInInverted, 543
- IsValveDigitalInInvertedEvent, 554
- IsValveOpen, 543
- IsValveOpenEvent, 554
- IsValveOpenInAnalogMode, 543
- IsValveOpenInAnalogModeEvent, 554
- IsValveOpenInDigitalMode, 544
- IsValveOpenInDigitalModeEvent, 555
- LoadValveTable, 544
- OnGetActiveRunningTableNumber, 544
- OnGetAnalogThresholdHigh, 544
- OnGetAnalogThresholdLow, 544
- OnGetAnalogVoltage, 544
- OnGetCurrentNumberOfValves, 545
- OnGetDigitalOutPortValve, 545
- OnGetDigitalPortDirection, 545
- OnGetDisplayMode, 545
- OnGetTableNamebyIndex, 545
- OnGetValveActive, 545
- OnGetValveBoardRevision, 545
- OnGetValveDigitalInPort, 545
- OnGetValveManualGroup, 545
- OnGetValveManualState, 546
- OnGetValveMode, 546
- OnIsDigitalOutPortInverted, 546
- OnIsValveDigitalInInverted, 546
- OnIsValveOpen, 546
- OnIsValveOpenInAnalogMode, 546
- OnIsValveOpenInDigitalMode, 546
- OnTableEntryChanged, 546
- SetActiveRunningTableNumber, 546
- SetAnalogThresholdHigh, 547
- SetAnalogThresholdLow, 547
- SetCurrentEditTableNumber, 547
- SetDefault, 547
- SetDigitalOutPortInvert, 547
- SetDigitalOutPortValve, 548
- SetDigitalPortDirection, 548
- SetDisplayMode, 548
- SetTableName, 549
- SetTableStep, 549
- SetTableStepAll, 549
- SetValveActive, 549
- SetValveDigitalInInvert, 550
- SetValveDigitalInPort, 550
- SetValveManualGroup, 550
- SetValveManualState, 550
- SetValveMode, 551
- SetValvesActiveMap, 551
- SetValvesManualStateMap, 551
- SetValveTableEntry, 551
- StoreValveTable, 552
- TableEntryChangedEvent, 555
- CWarnerValveControllerDeviceTesterFunctionNet, 555
 - !CWarnerValveControllerDeviceTesterFunctionNet, 556
 - ~CWarnerValveControllerDeviceTesterFunctionNet, 556
 - CWarnerValveControllerDeviceTesterFunctionNet, 556
 - GetIO, 556
 - GetSync, 556
 - SetADC, 557
 - SetIO, 557
 - SetIODirection, 557
 - SetTrigger, 557
 - SetTriggerSyncDirection, 558
- CWClassicFunctionNet, 558
 - CWClassicFunctionNet, 559
 - GetFilterParametersHeadstage, 559
 - GetHasChecksum, 559
 - GetHasRedLedHeadstage, 559
 - GetHeadstageOnOff, 559
 - GetResetFilter, 560
 - GetRFCConnectionStatus, 560

- GetRFFrequencyHeadstage, [560](#)
- GetRFFrequencyReceiver, [560](#)
- GetRFPower, [560](#)
- GetScanHeadstagesResult, [560](#)
- GetSelectedHeadstage, [560](#)
- GetSerialNumberHeadstage, [560](#)
- GetWPADebugMode, [560](#)
- GetWPAType, [560](#)
- ResetChannelmap, [561](#)
- ScanForHeadstages, [561](#)
- SetChannelmap, [561](#)
- SetFilterParametersHeadstage, [561](#)
- SetHasChecksum, [561](#)
- SetHeadstageOnOff, [561](#)
- SetHWSelectedChannels, [561](#)
- SetResetFilter, [561](#)
- SetRFFrequencyHeadstage, [561](#)
- SetRFFrequencyReceiver, [562](#)
- SetRFFrequencyReceiverEeprom, [562](#)
- SetRFLostBehaviour, [562](#)
- SetRFPower, [562](#)
- SetSelectedHeadstage, [562](#)
- SetSerialNumberHeadstage, [562](#)
- SetWPADebugMode, [562](#)
- SetWPAType, [562](#)
- CWirelessBaseFunctionNet, [563](#)
 - CreateWirelessHeadstageSerialNumberString, [563](#)
 - CWirelessBaseFunctionNet, [563](#)
- DacqDevice
 - CSafeISDeviceNet, [395](#)
- dacqgroup
 - CW2100_FunctionNet::AudioChannelsNet, [27](#)
- DACResolution
 - W2100_StimulusParametersNet, [594](#)
- DataState
 - HeadStageIDTypeState, [579](#)
- DeepCopy
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [343](#)
- DefineAmplification
 - CPgaDeviceNet, [315](#)
- DefineFrequencyRange
 - CPgaDeviceNet, [315](#)
- DefineNumAmplifications
 - CPgaDeviceNet, [315](#)
- DefineNumFrequencyRanges
 - CPgaDeviceNet, [315](#)
- DeselectAllHeadstages
 - CW2100_FunctionNet, [502](#)
- DeselectHeadstage
 - CW2100_FunctionNet, [502](#)
- DetectChipType
 - CCMOSMea_FunctionNet, [34](#)
- DEVICE_NOT_FOUND
 - Mcs::Usb, [25](#)
- DeviceArrival
 - CMcsUsbListNet, [223](#)
- DeviceData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- DeviceDataLength
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- DeviceId
 - CMcsUsbListEntryNet, [220](#)
- DeviceIdNet, [563](#)
 - BcdDevice, [565](#)
 - BusType, [565](#)
 - DeviceIdNet, [564](#)
 - IdProduct, [565](#)
 - IdVendor, [565](#)
 - operator=, [564](#)
- DeviceName
 - CMcsUsbListEntryNet, [220](#)
- DeviceRemoval
 - CMcsUsbListNet, [223](#)
- DigitalSource
 - DigitalSource< digitalsourceenum >, [565](#)
- DigitalSource< digitalsourceenum >, [565](#)
 - DigitalSource, [565](#)
 - MaxBitNumber, [565](#)
 - size, [566](#)
 - Source, [566](#)
- DisableAutoReset
 - CStg200xDownloadBasicNet, [446](#)
- DisableMultiFileMode
 - CStg200xDownloadNet, [455](#)
- Disconnect
 - CMcsUsbNet, [230](#)
- DisConnectDevice
 - CRadioControlledDevicesNet, [342](#)
- DoRamp
 - CRoboDacqNet, [357](#)
- DownloadFirmware
 - CMcsUsbFactoryNet, [209](#)
- DriverVersionNet, [566](#)
 - ~DriverVersionNet, [567](#)
 - DriverVersionNet, [567](#)
 - DriverVersionNet::FormatVersion, [567](#)
 - GetDestinationCode, [567](#)
 - GetDestinationName, [567](#), [569](#)
 - GetMajor, [569](#)
 - GetMinor, [569](#), [570](#)
 - GetNumEntries, [570](#)
 - GetStatus, [570](#)
 - GetVersionInt, [570](#), [571](#)
 - GetVersionString, [571](#)
- DriverVersionNet::FormatVersion
 - DriverVersionNet, [567](#)
- DSP
 - FirmwareDestinationNames, [572](#)
- DummyCommand
 - CLIH3DeviceNet, [109](#)
- Duration
 - CStimulusFunctionNet::SidebandData, [591](#)

- ElectricalStimulation
 - HeadStageIDType, [576](#)
- EmptyKey
 - CMcsUsbNet, [230](#)
- Emu_GetCellCapacity
 - CRoboDacqNet, [357](#)
- Emu_GetCellPotential
 - CRoboDacqNet, [357](#)
- Emu_GetCellResists
 - CRoboDacqNet, [357](#)
- Emu_GetElectrodeResists
 - CRoboDacqNet, [357](#)
- Emu_GetNoise
 - CRoboDacqNet, [357](#)
- Emu_SetCellCapacity
 - CRoboDacqNet, [357](#)
- Emu_SetCellPotential
 - CRoboDacqNet, [357](#)
- Emu_SetCellResists
 - CRoboDacqNet, [357](#)
- Emu_SetElectrodeResists
 - CRoboDacqNet, [358](#)
- Emu_SetNoise
 - CRoboDacqNet, [358](#)
- EnableAnalogOut
 - CSCUFunctionNet, [400](#)
- EnableAutoReset
 - CStg200xDownloadBasicNet, [447](#)
- EnableChannelsInGroup
 - CCMOSMea_FunctionNet, [34, 35](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [51, 52](#)
- EnableChecksum
 - CMeaDeviceNet, [262](#)
 - COctoPotDeviceNet, [304](#)
- EnableDigitalIn
 - CMeaDeviceNet, [262, 263](#)
 - COctoPotDeviceNet, [304](#)
- EnableExceptions
 - CMcsUsbNet, [230](#)
- EnableMultiFileMode
 - CStg200xDownloadNet, [455](#)
- EnableTimestamp
 - CMeaDeviceNet, [263](#)
 - COctoPotDeviceNet, [304](#)
- EnableUserTrigger
 - CLIH3DeviceNet, [110](#)
- enCMosMeaChipType
 - Mcs::Usb, [25](#)
- EnSTG200x_STATUS
 - Mcs::Usb, [25](#)
- Entry
 - HeadStageIDType, [576](#)
- Equals
 - CFilterCoefficientsNet, [65](#)
 - CMcsUsbListEntryNet, [218](#)
- HeadStageIDType, [576](#)
- HeadstageIDTypeObject, [578](#)
- EraseEepromRegisterPreconfig
 - CMcsUsbNet, [230](#)
- EraseFilterParameterPermanent
 - CFilterConfigurationNet, [67](#)
 - CFilterConfigurationRegisterNet, [69](#)
- Error_Callback_Aquisition_Stopped
 - CMcsUsbDacqNet, [204](#)
- Error_Callback_Data_lost
 - CMcsUsbDacqNet, [204](#)
- Error_Callback_Frames_Lost
 - CMcsUsbDacqNet, [204](#)
- Error_Callback_Packet_Error
 - CMcsUsbDacqNet, [204](#)
- Error_Callback_Queue_Full
 - CMcsUsbDacqNet, [204](#)
- Error_Callback_RingQueue_Full
 - CMcsUsbDacqNet, [204](#)
- ErrorEvent
 - CMcsUsbDacqNet, [205](#)
- FactoryReset
 - CTcxDeviceNet, [476](#)
- FeedbackGetSampleTimerCount
 - CMeaFeedbackFunctionNet, [271](#)
- FeedbackSetAnalogSource
 - CMeaFeedbackFunctionNet, [271](#)
- FeedbackSetChannelFilter
 - CMeaFeedbackFunctionNet, [271](#)
- FeedbackSetDigitalMapping
 - CMeaFeedbackFunctionNet, [271](#)
- FeedbackSetFeedback
 - CMeaFeedbackFunctionNet, [271](#)
- FeedbackSetFilterOff
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetFilterParameter
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetFilterParameter32
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetGlobalChannelFilter
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetIIRFilterParameter
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetLogic
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetMkFilter
 - CMeaFeedbackFunctionNet, [272](#)
- FeedbackSetNumberOfLogics
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetNumberOfRateCounter
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetNumberOfRateDetectors
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetNumberOfSpikeDetectors
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetNumberOfTriggers
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetRateCounter

- CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetRateDetector
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetSpikeDetectorThreshold
 - CMeaFeedbackFunctionNet, [273](#)
- FeedbackSetTrigger
 - CMeaFeedbackFunctionNet, [274](#)
- FilterActive
 - CFilterPropertyNet, [71](#)
- FilterBand
 - CFilterPropertyNet, [71](#)
- FilterFamily
 - CFilterPropertyNet, [71](#)
- FilterType
 - CFilterPropertyNet, [71](#)
- FindFilter
 - CCreateFilterNet, [47](#)
- FindFirmwareVersionMagicInBuffer
 - CMcsUsbFactoryNet, [209](#)
- FindReference
 - CRoboDeviceNet, [371](#)
- FindReferencel
 - CRoboStatorDeviceNet, [388](#)
- FindReferencePhase0
 - CRoboDeviceNet::RoboMainLowLevelCommands, [584](#)
- FindReferencePhase0XY
 - CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands, [590](#)
- FindReferenceXY
 - CRoboStatorDeviceNet, [388](#), [389](#)
- FindReferenceZ
 - CRoboStatorDeviceNet, [389](#)
- FirePressurePulse
 - CPPCFunctionNet, [325](#)
- FirmwareDestinationNames, [571](#)
 - Altera, [572](#)
 - Bootstrap, [572](#)
 - BUS1_MCSBUS1, [572](#)
 - BUS1_MCSBUS2, [572](#)
 - DSP, [572](#)
 - FPGA2, [572](#)
 - FPGA3, [573](#)
 - FPGA4, [573](#)
 - FPGA5, [573](#)
 - FPGA6, [573](#)
 - MCSBUS1, [573](#)
 - MCSBUS10, [573](#)
 - MCSBUS11, [573](#)
 - MCSBUS12, [573](#)
 - MCSBUS13, [573](#)
 - MCSBUS2, [573](#)
 - MCSBUS3, [573](#)
 - MCSBUS4, [574](#)
 - MCSBUS5, [574](#)
 - MCSBUS6, [574](#)
 - MCSBUS7, [574](#)
 - MCSBUS8, [574](#)
 - MCSBUS9, [574](#)
 - MCU1, [574](#)
 - PIC, [574](#)
 - PIC2, [574](#)
 - PIC3, [574](#)
 - PIC4, [574](#)
 - USB, [575](#)
- FluidControlDevice
 - CSafeISDeviceNet, [396](#)
- ForceStatusEvent
 - CStg200xDownloadBasicNet, [447](#)
 - CStimulusFunctionNet, [463](#)
- FPGA2
 - FirmwareDestinationNames, [572](#)
- FPGA3
 - FirmwareDestinationNames, [573](#)
- FPGA4
 - FirmwareDestinationNames, [573](#)
- FPGA5
 - FirmwareDestinationNames, [573](#)
- FPGA6
 - FirmwareDestinationNames, [573](#)
- FromIntPtr
 - StgStatusNet, [592](#)
- FromPtr
 - StgStatusNet, [592](#)
- FX3MCSDDataAddress
 - CMcsUsbFactoryNet, [215](#)
- FX3MCSDDataDeviceIdOffset
 - CMcsUsbFactoryNet, [215](#)
- FX3MCSDDataVersionOffset
 - CMcsUsbFactoryNet, [215](#)
- FYIProgram
 - CFYIDeviceNet, [79](#)
- FYITemp
 - CFYIDeviceNet, [79](#)
- Gain
 - CMeaDeviceNet, [267](#)
- Get2AnalogInput
 - CMcsBus_SensorNet, [140](#)
- Get2DigitalInput
 - CMcsBus_SensorNet, [140](#)
- Get4ADC
 - CMcsBus_SensorNet, [140](#)
- Get4ADCAverage
 - CMcsBus_SensorNet, [140](#)
- Get4ADCCatchampAverageShift
 - CMcsBus_SensorNet, [140](#)
- Get4ADCMode
 - CMcsBus_SensorNet, [140](#)
- Get4DAC
 - CMcsBus_SensorNet, [140](#)
- GetAbsMaxCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [298](#)
- GetAccelGyroCurrentRate
 - CW2100_FunctionNet, [502](#)
- GetAccelGyroDesiredRate
 - CW2100_FunctionNet, [502](#)

GetAccelGyroEnabled
 CW2100_FunctionNet, 502

GetAccelRange
 CW2100_FunctionNet, 502

GetActiveRunningTableNumber
 CWarnerValveControllerDeviceNet, 536

GetActiveRunningTableNumberEvent
 CWarnerValveControllerDeviceNet, 552

GetAdapterCode
 CMealImpedanceDeviceNet, 275
 CTEERFunctionNet, 489

GetAdapterType
 CMcsUsbDacqNet, 180

GetAdc
 CFluidControlDeviceNet, 73

GetAdcDataFormat
 CMcsUsbDacqNet, 180

GetADCInputOffset
 CCMOSMea_FunctionNet, 35

GetADCOffset
 CLIH3DeviceNet, 110

GetAdcOffset
 COctoPotDeviceNet, 304

GetAdcOffsetU1
 CTEERFunctionNet, 490

GetAdcOffsetU2
 CTEERFunctionNet, 490

GetADCs
 CMcsBus_SensorNet, 140

GetADCsLoop
 CMcsBus_SensorNet, 141

GetAdcZero
 CMcsUsbDacqNet, 180

GetAirpressure
 CRoboDeviceNet, 371

GetAirpressureLimit
 CRoboDeviceNet, 371

GetAirValve
 CRoboDeviceNet, 371

GetAllDigout
 CRoboDacqNet, 358

GetAmplification
 CPgaDeviceNet, 315

GetAmplitude_nA
 CTEERFunctionNet, 490

GetAnalogGain
 CMeaDeviceNet, 264

GetAnalogOutADCRange
 CSCUFunctionNet, 400

GetAnalogOutChannel
 CW2100_FunctionNet, 502

GetAnalogOutChannels
 CSCUFunctionNet, 400

GetAnalogOutDACRange
 CSCUFunctionNet, 400

GetAnalogOutFilter
 CW2100_FunctionNet, 502

GetAnalogRanges
 CStg200xBasicNet, 418

GetAnalogResolution
 CStg200xBasicNet, 418

GetAnalogThresholdHigh
 CWarnerValveControllerDeviceNet, 537

GetAnalogThresholdHighEvent
 CWarnerValveControllerDeviceNet, 552

GetAnalogThresholdLow
 CWarnerValveControllerDeviceNet, 537

GetAnalogThresholdLowEvent
 CWarnerValveControllerDeviceNet, 552

GetAnalogValueUnit
 CMcsUsbDacqNet, 180

GetAnalogVoltage
 CPPCFunctionNet, 326
 CPPS_FunctionNet, 333
 CWarnerValveControllerDeviceNet, 537

GetAnalogVoltageEvent
 CWarnerValveControllerDeviceNet, 552

GetAnalogVoltageRange
 CPPCFunctionNet, 326

GetAnalogVoltages
 CPPS_FunctionNet, 333

GetArraySize
 CMealImpedanceDeviceNet, 275

GetAudioChannels
 CMeaAudioFunctionNet, 248, 249
 CW2100_FunctionNet, 502

GetAutocalibrationDisabled
 CStg200xBasicNet, 419

GetAvailableBaseSamplerates
 CCMOSMeaDeviceNet, 44

GetAvailableChambers
 CWarnerUssingFunctionNet, 518

GetAvailableDeviceList
 CRFFFunctionNet, 347

GetAvailableDeviceListEx
 CRFFFunctionNet, 347

GetAvailableHeadstages
 CSCUFunctionNet, 400
 CW2100_FunctionNet, 502

GetAvailableHeadstagesEvent
 CSCUFunctionNet, 412

GetAvailableMemory
 CStg200xBasicNet, 419
 CStimulusFunctionNet, 463

GetAvailableSampleRates
 CMcsUsbDacqNet::CHWInfo, 102

GetAvailableStateList
 CRFFFunctionNet, 348

GetAvailableStateListEx
 CRFFFunctionNet, 348

GetAvailableVoltageRangesInMicroVolt
 CMcsUsbDacqNet::CHWInfo, 102

GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt
 CMcsUsbDacqNet::CHWInfo, 102

GetAxisConfig

- CRoboDeviceNet::RoboMainLowLevelCommands, 584
- GetAxisLED
 - CRoboocyte2DeviceNet, 386
- GetAxisParametersSignedEeprom
 - CMcsBus_AxisParametersNet, 115
- GetAxisParametersUnsignedEeprom
 - CMcsBus_AxisParametersNet, 115
- GetBaseFrequency
 - CRFFunctionNet, 348
- GetBaseSamplerate
 - CCMOSMeaDeviceNet, 44
- GetBath
 - CCMOSMea_FunctionNet, 35
- GetBathMode
 - CCMOSMea_FunctionNet, 35
- GetBatteryState
 - CW2100_FunctionNet, 503
- GetBatteryVoltage
 - CMultiBatteryChargerDeviceNet, 283
- GetBiQuad
 - CCreateFilterNet, 47
- GetBiQuads
 - CCreateFilterNet, 47
- GetBlankingEnable
 - CStg200xBasicNet, 419, 420
- GetBoardTemp
 - CTcxDeviceNet, 477
- GetBoostAlwaysOnMode
 - CW2100_StimulatorFunctionNet, 509
- GetBoostPreTime
 - CW2100_StimulatorFunctionNet, 509
- GetBubbleState
 - CPPS_FunctionNet, 333
- GetBubbleStatus
 - CMcsBus_SensorNet, 141
- GetBuffer
 - CGenericDevelopDeviceNet, 86
- GetBusAddress
 - CMcsBusNet, 154
- GetBusAddressEeprom
 - CMcsBusNet, 154
- GetByteBuffer
 - CGenericDevelopDeviceNet, 87
- GetBytesAvailable
 - CSerialPortNet, 413
- GetBytesPerSample
 - CTEERFunctionNet, 490
- GetCalibration
 - CTcxDeviceNet, 477
- GetCalibrationDecp
 - CTcxDeviceNet, 477
- GetCalibrationMax
 - CTcxDeviceNet, 477
- GetCalibrationMin
 - CTcxDeviceNet, 477
- GetCapacityC
 - CRoboDacqNet, 358
- GetCapacityV
 - CRoboDacqNet, 358
- GetCapacityX
 - CRoboDacqNet, 358
- GetCardinalDacqSamplerate
 - CInterfaceboardFunctionNet, 107
- GetCardinalStgOutputrate
 - CInterfaceboardFunctionNet, 107
- GetChannel
 - CSw2to64DeviceNet, 472
- GetChannelDataFillSize
 - CMcsUsbDacqNet, 180
- GetChannelDataI16
 - CCMOSMeaDeviceNet, 44
- GetChannelDataI32
 - CCMOSMeaDeviceNet, 44
- GetChannelDataUI16
 - CCMOSMeaDeviceNet, 44
- GetChannelDataUI32
 - CCMOSMeaDeviceNet, 45
- GetChannelLayout
 - CMcsUsbDacqNet, 180
- GetChannels
 - CMultiBatteryChargerDeviceNet, 283
 - CSw2to64DeviceNet, 473
- GetChannelsCountOfChamber
 - CWarnerUssingFunctionNet, 518
- GetChannelsInBlock
 - CMcsUsbDacqNet, 181
- GetChannelState
 - CMultiBatteryChargerDeviceNet, 283
- GetChargeCapacity
 - CMultiBatteryChargerDeviceNet, 284
- GetChargeCurrent
 - CMultiBatteryChargerDeviceNet, 284
- GetChargingMode
 - CMultiBatteryChargerDeviceNet, 284
- GetChargingPCoefficient
 - CMultiBatteryChargerDeviceNet, 284
- GetChecksumFromFX3Image
 - CMcsUsbFactoryNet, 210
- GetCheckVoltage
 - CokuvisionStimulatorDeviceNet, 307
- GetClampAmpSerialNumber
 - CRoboDacqNet, 358
- GetClampMode
 - CTEERFunctionNet, 490
 - CWarnerUssingFunctionNet, 519
- GetCMOSDataDictionary
 - CCMOSMeaDeviceNet, 45
- GetCoilCommunication
 - CPositionIIDeviceNet, 317
- GetColorRgb
 - CMultiwellOptoStimFunctionNet, 298
- GetColorStr
 - CMultiwellOptoStimFunctionNet, 299
- GetCommand
 - CMcsBusNet, 154, 155

- CPedoterDeviceNet, 312
- CRoboDacqNet, 358
- GetComplianceVoltage
 - CWarnerUssingFunctionNet, 519
- GetConfiguration
 - CMcsUsbNet, 231
- GetConfigurationBit
 - CRoboDacqNet, 358
- GetConfigurationBitAxc
 - CRoboDacqNet, 358
- GetConfigurationBitBlu_Led
 - CRoboDacqNet, 359
- GetConfigurationBitBlu_LedToggleFast
 - CRoboDacqNet, 359
- GetConfigurationBitBlu_LedToggleSlow
 - CRoboDacqNet, 359
- GetConfigurationBitCC_Gen
 - CRoboDacqNet, 359
- GetConfigurationBitCV_Gen
 - CRoboDacqNet, 359
- GetConfigurationBitRC_Gen
 - CRoboDacqNet, 359
- GetConfigurationBitRed_Led
 - CRoboDacqNet, 359
- GetConfigurationBitRed_LedSaturation
 - CRoboDacqNet, 359
- GetConfigurationBitRed_LedToggleFast
 - CRoboDacqNet, 359
- GetConfigurationBitRed_LedToggleSlow
 - CRoboDacqNet, 359
- GetConfigurationBitRelais
 - CRoboDacqNet, 359
- GetConfigurationBitRV_Gen
 - CRoboDacqNet, 360
- GetConfigurationBits
 - CRoboDacqNet, 360
- GetConfigurationBitStream
 - CRoboDacqNet, 360
- GetConfigurationBitSupply
 - CRoboDacqNet, 360
- GetConnectedDevice
 - CRFFFunctionNet, 349
- GetControllerParams
 - CTEERFunctionNet, 490
- GetCrossTalkOffset
 - CRoboDacqNet, 360
- GetCrossTalkOptimum
 - CRoboDacqNet, 360
- GetCurrent
 - CTcxDeviceNet, 477
- GetCurrentAirvalve
 - CRoboDeviceNet, 371
- GetCurrentAirvalveLimit
 - CRoboDeviceNet, 372
- GetCurrentCycle
 - CMeaCoatDeviceNet, 256
- GetCurrentEditTableNumber
 - CWarnerValveControllerDeviceNet, 537
- GetCurrentEnable
 - CTEERFunctionNet, 491
- GetCurrentFactor
 - COkuvisionStimulatorDeviceNet, 307
- GetCurrentNumberOfValves
 - CWarnerValveControllerDeviceNet, 538
- GetCurrentNumberOfValvesEvent
 - CWarnerValveControllerDeviceNet, 552
- GetCurrentPosition
 - CRoboDeviceNet, 372
- GetCurrentPositionI
 - CRoboStatorDeviceNet, 389
- GetCurrentPositionXY
 - CRoboStatorDeviceNet, 389
- GetCurrentPositionZ
 - CRoboStatorDeviceNet, 389
- GetCurrentRangeInNanoAmp
 - CStg200xBasicNet, 420
 - CStimulusFunctionNet, 464
 - CW2100_StimulatorFunctionNet, 509
- GetCurrentResolutionInNanoAmp
 - CStg200xBasicNet, 420
 - CStimulusFunctionNet, 464
 - CW2100_StimulatorFunctionNet, 509
- GetCycle
 - CMeaCleanDeviceNet, 252
- GetCycles
 - CMeaCleanDeviceNet, 252
 - CMeaCoatDeviceNet, 256
- GetD
 - CTcxDeviceNet, 477
- GetDacAmplificationFactor
 - CStg200xBasicNet, 421
- GetDacIdleValue
 - CLIH3DeviceNet, 110
- GetDacNampsPerDigitHighCurrentRange
 - CWarnerUssingFunctionNet, 519
- GetDACOffset
 - COkuvisionStimulatorDeviceNet, 308
- GetDacOffset
 - CDacCalibrationFunctionNet, 49
 - COctoPotDeviceNet, 304
- GetDacPampsPerDigitLowCurrentRange
 - CWarnerUssingFunctionNet, 520
- GetDacqRunStatus
 - CLIH3DeviceNet, 110
- GetDacRange
 - CW2100_FunctionNet, 503
- GetDACResolution
 - CStg200xBasicNet, 421
 - CStimulusFunctionNet, 464
 - CW2100_StimulatorFunctionNet, 509
- GetDACs
 - CMcsBus_SensorNet, 141
- GetDacUseldleValue
 - CLIH3DeviceNet, 111
- GetDacZero
 - CTEERFunctionNet, 491

- CWarnerUssingFunctionNet, 520
- GetDataFormat
 - CMcsUsbDacqNet, 181
- GetDataMode
 - CMcsUsbDacqNet, 181
- GetDDecp
 - CTcxDeviceNet, 477
- GetDestination
 - CMcsUsbFactoryNet, 210
- GetDestinationCode
 - DriverVersionNet, 567
- GetDestinationDisplayLabel
 - CMcsUsbFactoryNet, 210
- GetDestinationName
 - CMcsUsbFactoryNet, 210
 - DriverVersionNet, 567, 569
- GetDestinationSerialNumber
 - CMcsUsbFactoryNet, 210
- GetDestinationTargetAddress
 - CMcsUsbFactoryNet, 210
- GetDetectionThreshold
 - CMcsBus_SensorNet, 141
- GetDetectorValue
 - CMcsBus_SensorNet, 141
- GetDevice
 - CTcxDeviceNet, 477
- GetDeviceCannotStallOutRequests
 - CMcsUsbNet, 231
- GetDeviceCapableSpeed
 - CMcsUsbNet, 231
- GetDeviceEnum
 - CMcsUsbNet, 231
- GetDeviceGroupChannelInfos
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, 52
- GetDeviceId
 - CMcsUsbNet, 231
- GetDeviceList
 - CPositionImpDeviceNet, 321
- GetDeviceNames
 - CRadioControlledDevicesNet, 342
- GetDeviceRootHubVendorEnum
 - CMcsUsbNet, 231
- GetDeviceRootHubVendorID
 - CMcsUsbNet, 231
- GetDeviceRootHubVendorName
 - CMcsUsbNet, 231
- GetDeviceSpeed
 - CMcsUsbNet, 232
- GetDeviceType
 - CTcxDeviceNet, 478
- GetDevname
 - CTcxDeviceNet, 478
- GetDigin
 - CFluidControlDeviceNet, 74
- GetDiginState
 - CLIH3DeviceNet, 111
- GetDiginValue
 - CStg200xBasicNet, 421
- GetDigitalData
 - CMeaDigitalDataFunctionNet, 269
- GetDigitalIn
 - CPPCFunctionNet, 326
 - CPPS_FunctionNet, 333
- GetDigitalOutPortValve
 - CWarnerValveControllerDeviceNet, 538
- GetDigitalOutPortValveEvent
 - CWarnerValveControllerDeviceNet, 553
- GetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, 538
- GetDigitalPortDirectionEvent
 - CWarnerValveControllerDeviceNet, 553
- GetDigitalSource
 - CMcsUsbDacqNet, 181–183
- GetDigitalStimulatorTrigger
 - CW2100_StimulatorFunctionNet, 509
- GetDigitalStimulatorTriggerSlope
 - CW2100_StimulatorFunctionNet, 510
- GetDigout
 - CFluidControlDeviceNet, 74
 - CRoboDacqNet, 360
- GetDigoutMode
 - CStg200xBasicNet, 421
- GetDigoutValue
 - CStg200xBasicNet, 422
- GetDIO
 - CMcsBus_FYIExtensionNet, 118
- GetDischargeCapacity
 - CMultiBatteryChargerDeviceNet, 285
- GetDischargeCurrent
 - CMultiBatteryChargerDeviceNet, 285
- GetDischargeCurrentSetPoint
 - CMultiBatteryChargerDeviceNet, 285
- GetDisplayMode
 - CWarnerValveControllerDeviceNet, 538
- GetDisplayModeEvent
 - CWarnerValveControllerDeviceNet, 553
- GetDisplayText
 - CRoboDacqNet, 360
- GetDMax
 - CTcxDeviceNet, 478
- GetDMin
 - CTcxDeviceNet, 478
- GetDownsampleFactor
 - CRoboDacqNet, 360
- GetDSPHighPassByIndex
 - CIntanMea_FunctionNet, 104
- GetDuration
 - CMeaCoatDeviceNet, 256
- GetDuty
 - CTcxDeviceNet, 478
- GetEEPromPage
 - CLIH3DeviceNet, 111
- GetElectrodeDacMux

- CStg200xBasicNet, [422](#)
- GetElectrodeEnable
 - CStg200xBasicNet, [423](#)
- GetElectrodeMode
 - CStg200xBasicNet, [424](#)
- GetEnableAmplifierProtectionSwitch
 - CStg200xBasicNet, [424](#), [425](#)
- GetEnabledChannelsInGroup
 - CCMOSMea_FunctionNet, [35](#)
 - CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >, [52](#)
- GetEnableHeaterLimit
 - CTcxDeviceNet, [478](#)
- GetEnableThermocouple
 - CTcxDeviceNet, [478](#)
- GetEntry
 - CMcsUsbListEntryNet, [218](#), [219](#)
- GetEntryCount
 - CMcsUsbListEntryNet, [219](#)
- GetEnumerationSpeed
 - CMeaDeviceNet, [264](#)
- GetErrorAirpressure
 - CRoboDeviceNet, [372](#)
- GetErrorCurrentAirvalve
 - CRoboDeviceNet, [372](#)
- GetErrorText
 - CMcsUsbNet, [232](#)
- GetErrorVoltage12V
 - CRoboDeviceNet, [372](#)
- GetErrorVoltage5V
 - CRoboDeviceNet, [372](#)
- GetErrorVoltageAirvalve
 - CRoboDeviceNet, [372](#)
- GetErrorVoltageRs485A
 - CRoboDeviceNet, [372](#)
- GetErrorVoltageRs485B
 - CRoboDeviceNet, [373](#)
- GetErrorVoltageValves
 - CRoboDeviceNet, [373](#)
- GetExternalElectrodeEnable
 - CStg200xBasicNet, [425](#)
- GetFAAmplification
 - CStg200xBasicNet, [426](#)
- GetFilter
 - CRoboDacqNet, [360](#)
- GetFilterCoeffs
 - CRoboDacqNet, [360](#)
- GetFilterParametersHeadstage
 - CWClassicFunctionNet, [559](#)
- GetFilterProperties
 - CSCUFunctionNet, [401](#)
 - CW2100_FunctionNet, [503](#)
- GetFilterProperty
 - CMcsUsbDacqNet, [183](#)
 - CSCUFunctionNet, [401](#)
 - CW2100_FunctionNet, [503](#)
- GetFinalDischargeVoltage
 - CMultiBatteryChargerDeviceNet, [286](#)
- GetFirmwareVersion
 - CMcsUsbNet, [232](#)
- GetFirmwareVersionFromFile
 - CMcsUsbFactoryNet, [211](#)
- GetFirmwareVersionFromHexFile
 - CMcsUsbFactoryNet, [211](#)
- GetFPGA FirmwareType
 - CW2100_FunctionNet, [503](#)
- GetFrequency
 - CRadioControlledDevicesNet, [342](#)
- GetFrequencyRange
 - CPgaDeviceNet, [316](#)
- GetGain
 - CMeaDeviceNet, [264](#)
 - CPgaDeviceNet, [316](#)
- GetGate
 - CCMOSMea_FunctionNet, [35](#)
- GetGilsonDevice
 - CRoboocyte2DeviceNet, [386](#)
- GetGlobalRepeat
 - CDigOutStimulatorFunctionNet, [58](#)
- GetGNDI
 - CCMOSMea_FunctionNet, [35](#)
- GetGroupADCBits
 - CCMOSMea_FunctionNet, [35](#)
- GetGroupChannelBitmaskBySelect
 - CCMOSMea_FunctionNet, [36](#)
- GetGroupChannelBitmaskHS1NCBathCurrent
 - CCMOSMea_FunctionNet, [36](#)
- GetGroupChannelBitmaskHS1NCCol2Current
 - CCMOSMea_FunctionNet, [36](#)
- GetGroupChannelBitmaskHS1NChipTemp
 - CCMOSMea_FunctionNet, [36](#)
- GetGroupChannelBitmaskHS1Sidebands
 - CCMOSMea_FunctionNet, [37](#)
- GetGroupChannelBitmaskHS1TriggerStatus
 - CCMOSMea_FunctionNet, [37](#)
- GetGroupChannelBitmaskIFDigChannels
 - CCMOSMea_FunctionNet, [37](#)
- GetGroupChannelBitmaskInterfaceADC
 - CCMOSMea_FunctionNet, [37](#)
- GetGroupChannelBitmaskPacketFrameContext
 - CCMOSMea_FunctionNet, [38](#)
- GetGroupChannelBitmaskSTG1DACSignal
 - CCMOSMea_FunctionNet, [38](#)
- GetGroupChannelDataI16
 - CMcsUsbDacqNet, [183](#)
- GetGroupChannelDataI32
 - CMcsUsbDacqNet, [184](#)
- GetGroupChannelDataUI16
 - CMcsUsbDacqNet, [184](#)
- GetGroupChannelDataUI32
 - CMcsUsbDacqNet, [185](#)
- GetGroupDCOffset
 - CCMOSMea_FunctionNet, [38](#)
- GetGroupID

- CCMOSMea_FunctionNet, [38](#)
- CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [52](#)
- GetGroupNumberOfChannels
 - CCMOSMea_FunctionNet, [39](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [52](#)
- GetGroupResolutionPerDigit
 - CCMOSMea_FunctionNet, [39](#)
- GetGroupSampleSize
 - CCMOSMea_FunctionNet, [39](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [53](#)
- GetGroupType
 - CCMOSMea_FunctionNet, [39](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [53](#)
- GetGroupUnit
 - CCMOSMea_FunctionNet, [39](#), [40](#)
- GetGyroRange
 - CW2100_FunctionNet, [503](#)
- GetHardwareMaxRange
 - CMcsUsbDacqNet, [186](#)
- GetHardwareMinRange
 - CMcsUsbDacqNet, [186](#)
- GetHardwareRevision
 - CMcsUsbNet, [232](#)
- GetHasChecksum
 - CWClassicFunctionNet, [559](#)
- GetHashCode
 - HeadstageIDTypeObject, [578](#)
- GetHasRedLedHeadstage
 - CWClassicFunctionNet, [559](#)
- GetHasThermocouple
 - CTcxDeviceNet, [478](#)
- GetHeadstage
 - CStg200xBasicNet, [426](#)
- GetHeadstageActive
 - CMcsUsbNet, [233](#)
- GetHeadstageAdcBits
 - CSCUFunctionNet, [401](#)
- GetHeadstageAdcRangeInMicroVolt
 - CSCUFunctionNet, [402](#)
- GetHeadstageDacBits
 - CSCUFunctionNet, [402](#)
- GetHeadstageDacCurrentRangeInMicroAmpere
 - CSCUFunctionNet, [402](#)
- GetHeadstageDacCurrentResolutionInNanoAmpere
 - CSCUFunctionNet, [403](#)
- GetHeadstageDacVoltageRangeInMilliVolt
 - CSCUFunctionNet, [403](#)
- GetHeadstageDacVoltageResolutionInMicroVolt
 - CSCUFunctionNet, [403](#)
- GetHeadstageGainInPer mille
 - CSCUFunctionNet, [404](#)
- GetHeadstageID
 - CMcsUsbNet, [233](#)
 - CSCUFunctionNet, [404](#)
- GetHeadstageNumberOfAnalogChannels
 - CSCUFunctionNet, [404](#)
- GetHeadstageNumberOfStimulationChannels
 - CSCUFunctionNet, [405](#)
- GetHeadstageOnOff
 - CW2100_FunctionNet, [503](#)
 - CWClassicFunctionNet, [559](#)
- GetHeadstagePresent
 - CMcsUsbNet, [233](#)
- GetHeadstageSamplerate
 - CSCUFunctionNet, [405](#)
- GetHeadstageSamplingActive
 - CW2100_FunctionNet, [503](#)
- GetHeadstageSerialNumber
 - CSCUFunctionNet, [405](#)
- GetHeaterLimit
 - CTcxDeviceNet, [479](#)
- GetHeaterTemp
 - CTcxDeviceNet, [479](#)
- GetHighCurrentRange
 - CWarnerUssingFunctionNet, [520](#)
- GetHighpassFilterEnable
 - CFilterConfigurationNet, [67](#)
- GetHWConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands,
[584](#)
- GetHWRevision
 - CRoboDeviceNet::RoboMainLowLevelCommands,
[584](#)
- GetHWRevisionEeprom
 - CMcsBusNet, [155](#)
- GetI
 - CTcxDeviceNet, [479](#)
- GetIC
 - CRoboDacqNet, [361](#)
- GetIClamp
 - CRoboDacqNet, [361](#)
- GetICoeff
 - CRobo_FYITemp_FunctionNet, [353](#)
- GetICOffset
 - CRoboDacqNet, [361](#)
- GetIDecp
 - CTcxDeviceNet, [479](#)
- GetIDent
 - CMcsUsbNet, [233](#)
- GetIdleModeOffset
 - CWarnerUssingFunctionNet, [521](#)
- GetIGain
 - CRoboDacqNet, [361](#)
- GetIMax

- CTcxDeviceNet, 479
- GetIMin
 - CTcxDeviceNet, 479
- GetImpedanceResult
 - CIntanMea_FunctionNet, 104
- GetImpedanceTestFrequency
 - CMealImpedanceDeviceNet, 275
- GetImpId
 - CPositionImpDeviceNet, 321
- GetImplantCurrentSetpoint
 - CPositionIIDeviceNet, 318
- GetImplantState
 - CPositionIIDeviceNet, 318
- GetInMovement
 - CRoboDeviceNet, 373
- GetIntanRegister
 - CIntanMea_FunctionNet, 104
- GetIntBuffer
 - CGenericDevelopDeviceNet, 87
- GetIO
 - CWarnerValveControllerDeviceTesterFunctionNet, 556
- GetIOut
 - CTcxDeviceNet, 479
- GetLastAnswer
 - CGilsonDeviceNet, 98
- GetLastUSBError
 - CMcsUsbNet, 234
- GetLatency
 - CMcsBus_SensorNet, 141
- GetLatencyCounter
 - CMcsBus_SensorNet, 141
- GetLayoutConfiguration
 - CMEA2100x256FunctionNet, 246
- GetLEDSwitch
 - CMcsBus_ExtensionNet, 117
- GetLength
 - CRobo_FYIProgram_FunctionNet, 351
- GetLiquidResistance
 - CTEERFunctionNet, 491
 - CWarnerUssingFunctionNet, 521
- GetListModelIndexRange
 - CStg200xBasicNet, 426
- GetListmodeTriggerSource
 - CStg200xBasicNet, 426
- GetLowCurrentRange
 - CWarnerUssingFunctionNet, 521
- GetLowerFrequencyByIndex
 - CIntanMea_FunctionNet, 104
- GetMajor
 - DriverVersionNet, 569
- GetMaxChunkSize_Byte
 - CTEERFunctionNet, 491
- GetMaxCurrent
 - CMeaCoatDeviceNet, 256
- GetMaxDurationHighCurrentInMicroSec
 - CMultiwellOptoStimFunctionNet, 299
- GetMaxDutyCycleHighCurrent
 - CMultiwellOptoStimFunctionNet, 299
- GetMaxHeaterPowerMultiwell
 - CTcxDeviceNet, 480
- GetMaxNoPressure
 - CRoboDeviceNet::RoboMainLowLevelCommands, 584
- GetMaxNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, 584
- GetMaxNumberOfHeadstages
 - CSCUFunctionNet, 406
- GetMaxNumOfColumns
 - CCMOSMea_FunctionNet, 40
- GetMaxP
 - CTcxDeviceNet, 480
- GetMaxpDecp
 - CTcxDeviceNet, 480
- GetMaxpMax
 - CTcxDeviceNet, 480
- GetMaxpMin
 - CTcxDeviceNet, 480
- GetMaxPower
 - COKuvisionStimulatorDeviceNet, 308
 - CRobo_FYITemp_FunctionNet, 353
- GetMaxPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, 584
- GetMaxSamplingFrequency
 - CMcsUsbDacqNet, 186
- GetMaxStimulusChannelsPerHeadstage
 - CSCUFunctionNet, 406
- GetMaxVoltage
 - CMeaCleanDeviceNet, 252
 - COKuvisionStimulatorDeviceNet, 308
- GetMCAcceleration
 - CMcsBus_MotorControlNet, 122
- GetMCAccelerationEeprom
 - CMcsBus_MotorControlNet, 122
- GetMCAccelerationShortCommand
 - CMcsBus_MotorControlNet, 122
- GetMCAXisRevisionEeprom
 - CMcsBus_MotorControlNet, 122
- GetMCBreakCurrent
 - CMcsBus_MotorControlNet, 122
- GetMCBreakCurrentEeprom
 - CMcsBus_MotorControlNet, 123
- GetMCCConfig
 - CMcsBus_MotorControlNet, 123
- GetMCCConfigEeprom
 - CMcsBus_MotorControlNet, 123
- GetMCCurrent
 - CMcsBus_MotorControlNet, 123
- GetMCCurrentEeprom
 - CMcsBus_MotorControlNet, 123
- GetMCCurrentMode
 - CMcsBus_MotorControlNet, 123
- GetMCCurrentModeEeprom
 - CMcsBus_MotorControlNet, 123

- GetMCCurrentModeShortCommand
 - CMcsBus_MotorControlNet, [124](#)
- GetMCCurrentPosition
 - CMcsBus_MotorControlNet, [124](#)
- GetMCCurrentShortCommand
 - CMcsBus_MotorControlNet, [124](#)
- GetMCCurrentSpeed
 - CMcsBus_MotorControlNet, [124](#)
- GetMCMaxAcceleration
 - CMcsBus_MotorControlNet, [124](#)
- GetMCMaxAccelerationEeprom
 - CMcsBus_MotorControlNet, [124](#)
- GetMCMaxCurrent
 - CMcsBus_MotorControlNet, [124](#)
- GetMCMaxCurrentEeprom
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMaxSpeed
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMaxSpeedEeprom
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMaxTravel
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMaxTravelEeprom
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMaxTravelShortCommand
 - CMcsBus_MotorControlNet, [125](#)
- GetMCMovement
 - CMcsBus_MotorControlNet, [125](#)
- GetMCNewPosition
 - CMcsBus_MotorControlNet, [126](#)
- GetMCOutputOnOff
 - CMcsBus_MotorControlNet, [126](#)
- GetMCPhase
 - CMcsBus_MotorControlNet, [126](#)
- GetMCPhaseOffset
 - CMcsBus_MotorControlNet, [126](#)
- GetMCReference
 - CMcsBus_MotorControlNet, [126](#)
- GetMCReferenceCurrent
 - CMcsBus_MotorControlNet, [126](#)
- GetMCReferenceCurrentEeprom
 - CMcsBus_MotorControlNet, [126](#)
- GetMCRegulatorGain
 - CMcsBus_MotorControlNet, [127](#)
- GetMCRegulatorGainEeprom
 - CMcsBus_MotorControlNet, [127](#)
- GetMcsBus_Extension
 - CRoboocyte2DeviceNet, [386](#)
- GetMCScalingFactor
 - CMcsBus_MotorControlNet, [127](#)
- GetMCScalingFactorEeprom
 - CMcsBus_MotorControlNet, [127](#)
- GetMCSpeed
 - CMcsBus_MotorControlNet, [127](#)
- GetMCSpeedEeprom
 - CMcsBus_MotorControlNet, [127](#)
- GetMCSpeedShortCommand
 - CMcsBus_MotorControlNet, [127](#)
- GetMCSpeedUnitEeprom
 - CMcsBus_MotorControlNet, [128](#)
- GetMCStandbyCurrent
 - CMcsBus_MotorControlNet, [128](#)
- GetMCStandbyCurrentEeprom
 - CMcsBus_MotorControlNet, [128](#)
- GetMCStandbyTime
 - CMcsBus_MotorControlNet, [128](#)
- GetMCStandbyTimeEeprom
 - CMcsBus_MotorControlNet, [128](#)
- GetMea21UsbPort
 - CMcsUsbNet, [234](#)
- GetMeaLayout
 - CMcsUsbDacqNet, [186](#)
- GetMemoryUsageDAC
 - CStg200xDownloadBasicNet, [447](#)
- GetMemoryUsageSyncout
 - CStg200xDownloadBasicNet, [447](#)
- GetMinimalThreshold
 - CMcsBus_SensorNet, [141](#)
- GetMinNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [584](#)
- GetMinor
 - DriverVersionNet, [569](#), [570](#)
- GetMinPressure
 - CRoboDeviceNet, [373](#)
 - CRoboDeviceNet::RoboMainLowLevelCommands, [584](#)
- GetMinPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [585](#)
- GetMinSamplingFrequencyStepsize
 - CMcsUsbDacqNet, [186](#)
- GetMinVoltage
 - CMeaCleanDeviceNet, [252](#)
- GetModeSelect
 - CPulseGeneratorFunctionNet, [339](#)
- GetModuleCurrent
 - CStg200xDownloadNet, [456](#)
- GetModuleTemp
 - CStg200xDownloadNet, [456](#)
- GetMovementError
 - CRoboDeviceNet, [373](#)
- GetMovePump
 - CMcsBus_SensorNet, [142](#)
- GetMultiHeadstageMode
 - CW2100_FunctionNet, [503](#)
- GetMultiplexedDataChannelsInBlock
 - CStimulusFunctionNet, [464](#)
- GetNanoVoltsPerKelvin
 - CMcsBus_TempSensorNet, [147](#)
- GetNeurochipMemoryData
 - CCMOSMea_FunctionNet, [40](#)
- GetNeurochipMemorySize
 - CCMOSMea_FunctionNet, [40](#)
- GetNIC_MS
 - CRoboDacqNet, [361](#)

- GetNUC_MS
 - CRoboDacqNet, [361](#)
- GetNumAmplifications
 - CPgaDeviceNet, [316](#)
- GetNumber
 - CMeaSwitchDeviceNet, [277](#)
 - CSw2to64DeviceNet, [473](#)
- GetNumberOfAnalogChannels
 - CStg200xBasicNet, [426](#)
 - CStimulusFunctionNet, [465](#)
 - CW2100_StimulatorFunctionNet, [510](#)
- GetNumberOfAudioChannels
 - CMeaAudioFunctionNet, [249](#)
- GetNumberOfAvailableChambers
 - CWarnerUssingFunctionNet, [522](#)
- GetNumberOfAvailableSamples
 - CTEERFunctionNet, [491](#)
- GetNumberOfChannels
 - CDigOutStimulatorFunctionNet, [58](#)
- GetNumberOfDataBits
 - CMcsUsbDacqNet, [187](#)
- GetNumberOfDevices
 - CMcsUsbListNet, [222](#)
- GetNumberOfHardwareSlotsForChambers
 - CWarnerUssingFunctionNet, [522](#)
- GetNumberOfHWADCCChannels
 - CMcsUsbDacqNet::CHWInfo, [102](#)
- GetNumberOfHWDACPaths
 - CStg200xBasicNet, [426](#)
- GetNumberOfHWDigitalChannels
 - CMcsUsbDacqNet::CHWInfo, [102](#)
- GetNumberOfStimulationElectrodes
 - CStg200xBasicNet, [426](#)
- GetNumberOfStimulationSourcesPerElectrode
 - CStg200xBasicNet, [427](#)
- GetNumberOfSupportedGroups
 - CCMOSMea_FunctionNet, [40](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-qGroupChannelEnumTemplateNet, Dac-qGroupChannelEnumTemplate, CDevice-GroupChannelInfoTemplateNet >, [53](#)
- GetNumberOfSyncoutChannels
 - CStg200xBasicNet, [427](#)
 - CW2100_StimulatorFunctionNet, [510](#)
- GetNumberOfTriggerInputs
 - CStg200xBasicNet, [427](#)
 - CW2100_StimulatorFunctionNet, [510](#)
- GetNumConfigurations
 - CMcsUsbNet, [234](#)
- GetNumControlChannels
 - CTcxDeviceNet, [480](#)
- GetNumDestinations
 - CMcsUsbFactoryNet, [211](#)
- GetNumDevices
 - CTcxDeviceNet, [480](#)
- GetNumEntries
 - DriverVersionNet, [570](#)
- GetNumFrequencyRanges
 - CPgaDeviceNet, [316](#)
- GetNumMeasureChannels
 - CTcxDeviceNet, [480](#)
- GetNUV_MS
 - CRoboDacqNet, [361](#)
- GetOffsetCurrent
 - CMeaCoatDeviceNet, [257](#)
- GetOnOff
 - CPositionIIDeviceNet, [318](#)
 - CTcxDeviceNet, [481](#)
- GetOutputCurrent
 - CMeaCoatDeviceNet, [257](#)
- GetOutputRate
 - CStg200xBasicNet, [427](#)
- GetOutputVoltage
 - CMeaCleanDeviceNet, [252](#)
- GetP
 - CTcxDeviceNet, [481](#)
- GetParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, [585](#)
- GetPattern
 - CMeaSwitchDeviceNet, [277](#)
- GetPatternBool
 - CMeaSwitchDeviceNet, [278](#)
- GetPauseDuration
 - CMeaCoatDeviceNet, [257](#)
- GetPCoeff
 - CRobo_FYITemp_FunctionNet, [353](#)
- GetPDecp
 - CTcxDeviceNet, [481](#)
- GetPeriod
 - CPulseGeneratorFunctionNet, [340](#)
- GetPeriod_us
 - CTEERFunctionNet, [492](#)
- GetPermanentCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [299](#)
- GetPGain
 - CRoboDacqNet, [361](#)
- GetPhases
 - CRoboDeviceNet::RoboMainLowLevelCommands, [585](#)
- GetPicFirmwareType
 - CW2100_FunctionNet, [503](#)
- GetPiezoState
 - CMcsBus_SensorNet, [142](#)
- GetPlateClampLockState
 - CMultiwellDeviceNet, [293](#)
- GetPlateClampState
 - CMultiwellDeviceNet, [293](#)
- GetPlateClampStateByHeadstage
 - CMultiwellCallbackFunctionNet, [289](#)
 - CMultiwellDeviceNet, [293](#)
- GetPlateClampStateByHeadstageEvent
 - CMultiwellCallbackFunctionNet, [291](#)
- GetPlateMux
 - CMultiwellDeviceNet, [294](#)
- GetPlateMuxByHeadstage

- CMultiwellDeviceNet, 294
- GetPlateType
 - CMultiwellDeviceNet, 294
- GetPlateTypeByHeadstage
 - CMultiwellCallbackFunctionNet, 290
 - CMultiwellDeviceNet, 294
- GetPlateTypeByHeadstageEvent
 - CMultiwellCallbackFunctionNet, 291
- GetPMax
 - CTcxDeviceNet, 481
- GetPMin
 - CTcxDeviceNet, 481
- GetPoti
 - CMcsUsbDacqNet, 187
- GetPOut
 - CTcxDeviceNet, 481
- GetPressure
 - CMcsBus_SensorNet, 142
- GetPressureOffset
 - CMcsBus_SensorNet, 142
- GetPressureRange
 - CPPCFunctionNet, 327
- GetPulseform
 - CokuvisionStimulatorDeviceNet, 308
- GetPulseLength
 - CPulseGeneratorFunctionNet, 340
- GetPumpCouple
 - CPPS_FunctionNet, 333
- GetPumpEnableSpeedRatio
 - CPPS_FunctionNet, 333
- GetPumpFastOnOff
 - CPPS_FunctionNet, 333
- GetPumpFastSpeed
 - CPPS_FunctionNet, 333
- GetPumpFunctionSpeeds
 - CPPS_FunctionNet, 333
- GetPumpManualOnOff
 - CPPS_FunctionNet, 333
- GetPumpMaxSpeed
 - CPPS_FunctionNet, 334
- GetPumpModeType
 - CPPCFunctionNet, 327
 - CPPS_FunctionNet, 334
- GetPumpSpeed
 - CRoboFluidDeviceNet, 383
- GetPumpSpeedRatio
 - CPPS_FunctionNet, 334
- GetPumpSpeedUnit
 - CPPCFunctionNet, 327
 - CPPS_FunctionNet, 334
- GetPWM
 - CFluidControlDeviceNet, 74
- GetPwrOut
 - CTcxDeviceNet, 481
- GetPwrSet
 - CTcxDeviceNet, 482
- GetRatedCapacity
 - CMultiBatteryChargerDeviceNet, 286
- GetReady
 - CMealImpedanceDeviceNet, 275
- GetRecordingNumber
 - CRoboDacqNet, 361
- GetReferenceElectrodeMode
 - CSCUFunctionNet, 406
- GetReferenceElectrodeSwitchState
 - CSCUFunctionNet, 406
- GetReferenceTemperature
 - CFluidControlDeviceNet, 74
- GetRegulationTimeouts
 - CMcsBus_SensorNet, 142
- GetRegulatorFactor
 - CMcsBus_SensorNet, 143
- GetRegulatorOnOff
 - CMcsBus_SensorNet, 143
 - CRobo_FYITemp_FunctionNet, 353
- GetRegulatorStatus
 - CMcsBus_SensorNet, 143
- GetRepeats
 - CProgramPressureCurveNet, 337
- GetRes1
 - CTcxDeviceNet, 482
- GetRes2
 - CTcxDeviceNet, 482
- GetResetFilter
 - CWClassicFunctionNet, 560
- GetResistanceC
 - CRoboDacqNet, 361
- GetResistanceV
 - CRoboDacqNet, 361
- GetResolutionPerDigit
 - CMcsUsbDacqNet, 187
- GetResS
 - CTcxDeviceNet, 482
- GetResult
 - CMealImpedanceDeviceNet, 275
- GetResX
 - CTcxDeviceNet, 482
- GetRFConnectionStatus
 - CWClassicFunctionNet, 560
- GetRFFrequency
 - CPositionImpDeviceNet, 321
- GetRFFrequencyHeadstage
 - CWClassicFunctionNet, 560
- GetRFFrequencyReceiver
 - CWClassicFunctionNet, 560
- GetRFPower
 - CWClassicFunctionNet, 560
- GetRoboDacq
 - CRoboocyte2DeviceNet, 386
- GetRoboFluidDevice
 - CEncapsulatorDeviceNet, 61
 - CRoboocyte2DeviceNet, 386
- GetRotaryPositionCode
 - CTEERFunctionNet, 492
- GetRotatePump
 - CMcsBus_SensorNet, 143

GetROut
 CTcxDeviceNet, [482](#)

GetRTC
 CokuvisionStimulatorDeviceNet, [308](#)

GetSampleBufferChunk
 CTEERFunctionNet, [492](#)

GetSampleInterval
 CLIH3DeviceNet, [112](#)

GetSamplePeriode
 CMcsBus_SensorNet, [143](#)

GetSamplerate
 CMcsUsbDacqNet, [187](#)

GetSampleVoltageBuffer_uV
 CTEERFunctionNet, [492](#)

GetScaleFactorU1
 CTEERFunctionNet, [493](#)

GetScaleFactorU2
 CTEERFunctionNet, [493](#)

GetScanHeadstagesResult
 CWClassicFunctionNet, [560](#)

GetScreen
 CRoboDacqNet, [362](#)

GetSearchReferenceFastAccel
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [585](#)

GetSearchReferenceFastSpeed
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [585](#)

GetSearchReferenceFineAccel
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [585](#)

GetSearchReferenceFineSpeed
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [585](#)

GetSearchReferenceMethod
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [586](#)

GetSearchReferenceMoveOut
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [586](#)

GetSearchReferenceOffsetPos
 CRoboDeviceNet::RoboMainLowLevelCommands,
 [586](#)

GetSelectedChannels
 CW2100_FunctionNet, [504](#)

GetSelectedHeadstage
 CWClassicFunctionNet, [560](#)

GetSelectedHeadstageState
 CW2100_FunctionNet, [504](#)

GetSensorType
 CTcxDeviceNet, [482](#)

GetSerialNumber
 CMcsUsbNet, [234](#)

GetSerialNumberHeadstage
 CWClassicFunctionNet, [560](#)

GetSetpoint
 CTcxDeviceNet, [482](#)

GetSetpointDecp
 CTcxDeviceNet, [483](#)

GetSetpointMax
 CTcxDeviceNet, [483](#)

GetSetpointMin
 CTcxDeviceNet, [483](#)

GetShortBuffer
 CGenericDevelopDeviceNet, [88](#)

GetSimulation
 CRoboDacqNet, [362](#)

GetSingleHeater
 CMcsBus_FYIExtensionNet, [118](#)

GetSingleValve
 CFluidControlDeviceNet, [75](#)
 CRoboFluidDeviceNet, [383](#)

GetSlope
 CMeaCleanDeviceNet, [253](#)
 CMeaCoatDeviceNet, [257](#)

GetSoftwareKey
 CMcsUsbNet, [234](#)

GetSoftwareKeyString
 CMcsUsbNet, [234](#)

GetSolIIPressure
 CMcsBus_SensorNet, [143](#)

GetSolITemp
 CRobo_FYITemp_FunctionNet, [353](#)

GetSourceBulk
 CCMOSMea_FunctionNet, [40](#)

GetSourceDrain
 CCMOSMea_FunctionNet, [40](#)

GetSourceGate
 CCMOSMea_FunctionNet, [40](#)

GetStartTriggerSlope
 CDigOutStimulatorFunctionNet, [58](#)

GetState
 CRFFFunctionNet, [349](#)
 CRobo_FYIProgram_FunctionNet, [351](#)

GetStatus
 CMcsUsbNet, [234](#)
 DriverVersionNet, [570](#)

GetStatusOfLastCommand
 CMcsUsbNet, [235](#)

GetStgProgramInfo
 CStg200xBasicNet, [427](#), [428](#)

GetStgVersionInfo
 CStg200xBasicNet, [428](#)

GetStimulationPatternMemory
 CW2100_StimulatorFunctionNet, [510](#)

GetStimulatorStatus
 CokuvisionStimulatorDeviceNet, [308](#)

GetStimulusParametersCache
 CW2100_FunctionNet, [504](#)

GetStimulusParametersFromSelectedHS
 CW2100_FunctionNet, [504](#)

GetStimulusSites
 CCMOSMea_FunctionNet, [41](#)

GetStiumlusParameters
 CW2100_FunctionNet, [504](#)

GetStopTriggerSlope

- CDigOutStimulatorFunctionNet, 59
- GetSubChannel
 - CMcsBus_MotorControlNet, 128
- GetSupplyVoltage
 - CPPCFunctionNet, 327
 - CPPS_FunctionNet, 334
- GetSweepCount
 - CStg200xDownloadBasicNet, 448
- GetSync
 - CWarnerValveControllerDeviceTesterFunctionNet, 556
- GetSyncoutMap
 - CStg200xBasicNet, 428
- GetSyncState
 - CMcsBus_SensorNet, 143
- GetTableName
 - CWarnerValveControllerDeviceNet, 539
- GetTableNamebyIndex
 - CWarnerValveControllerDeviceNet, 539
- GetTableNamebyIndexEvent
 - CWarnerValveControllerDeviceNet, 553
- GetTablepointer
 - CRetinalLedDeviceNet, 345
- GetTemperatur
 - CMcsBus_TempSensorNet, 147
- GetTestMode
 - CRFFFunctionNet, 349
- GetThermocoupleCalibration
 - CFluidControlDeviceNet, 75
 - CTcxDeviceNet, 483
- GetThermocoupleNanovoltPerKelvin
 - CFluidControlDeviceNet, 75
 - CTcxDeviceNet, 483
- GetThermocoupleReferenceTemp
 - CTcxDeviceNet, 483
- GetThermocoupleTemp
 - CTcxDeviceNet, 483
- GetThermocoupleTempAbs
 - CTcxDeviceNet, 484
- GetThermocoupleTemperature
 - CFluidControlDeviceNet, 76
- GetThermoOffset
 - CMcsBus_TempSensorNet, 147
- GetThermoTemp
 - CMcsBus_TempSensorNet, 148
- GetThermoVoltage
 - CMcsBus_TempSensorNet, 148
- GetTimeInPause
 - CMeaCoatDeviceNet, 257
- GetTimeInPlateau
 - CMeaCoatDeviceNet, 258
- GetTimeResolutionInNanoSeconds
 - CW2100_StimulatorFunctionNet, 510
- GetTimeSlot
 - CW2100_StimulatorFunctionNet, 510
- GetTotalMemory
 - CStg200xBasicNet, 429
 - CStimulusFunctionNet, 465
- GetTotalNumberOfDigitalPorts
 - CWarnerValveControllerDeviceNet, 539
- GetTotalNumberOfTables
 - CWarnerValveControllerDeviceNet, 539
- GetTotalNumberOfValves
 - CWarnerValveControllerDeviceNet, 539
- GetTotalTableSize
 - CWarnerValveControllerDeviceNet, 540
- GetTrigger
 - CStg200xDownloadBasicNet, 448
- GetTriggerSource
 - CStg200xBasicNet, 429
- GetU1Offset
 - CWarnerUssingFunctionNet, 522
- GetU1Reference
 - CWarnerUssingFunctionNet, 522
- GetU2Offset
 - CWarnerUssingFunctionNet, 523
- GetU2Reference
 - CWarnerUssingFunctionNet, 523
- GetUByteBuffer
 - CGenericDevelopDeviceNet, 89
- GetUC
 - CRoboDacqNet, 362
- GetUClamp
 - CRoboDacqNet, 362
- GetUCOffset
 - CRoboDacqNet, 362
- GetUIntBuffer
 - CGenericDevelopDeviceNet, 89
- GetUnit
 - CTcxDeviceNet, 484
- GetUnitDescription
 - CWarnerUssingFunctionNet, 523
- GetUnitExponent
 - CWarnerUssingFunctionNet, 524
- GetUnitName
 - CWarnerUssingFunctionNet, 524
- GetUnitsPerDigit
 - CWarnerUssingFunctionNet, 524
- GetUOut
 - CTcxDeviceNet, 484
- GetUpdateDisplay
 - CRoboDacqNet, 362
- GetUpperFrequencyByIndex
 - CIntanMea_FunctionNet, 105
- GetUptimeSeconds
 - CTEERFunctionNet, 493
 - CWarnerUssingFunctionNet, 526
- GetUSBDeviceIDFromFX3Image
 - CMcsUsbFactoryNet, 211
- GetUsbListEntries
 - CMcsUsbListNet, 222
- GetUsbListEntry
 - CMcsUsbListNet, 222
 - CMcsUsbNet, 235
- GetUseBubble
 - CPPS_FunctionNet, 334

- GetUserCodeFromBitFile
 - CMcsUsbFactoryNet, [211](#)
- GetUserCodeFromFlash
 - CMcsUsbFactoryNet, [211](#)
- GetUserDefinedName
 - CW2100_FunctionNet, [504](#)
- GetUserDefinedNameCache
 - CW2100_FunctionNet, [504](#)
- GetUserDefinedNameFromSelectedHS
 - CW2100_FunctionNet, [505](#)
- GetUserParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, [586](#)
- GetUShortBuffer
 - CGenericDevelopDeviceNet, [90](#)
- GetUV
 - CRoboDacqNet, [362](#)
- GetUVOffset
 - CRoboDacqNet, [362](#)
- GetValue
 - CTcxDeviceNet, [484](#)
- GetValueHires
 - CTcxDeviceNet, [484](#)
- GetValve
 - CFluidControlDeviceNet, [76](#)
 - CRoboFluidDeviceNet, [383](#)
- GetValve1
 - CRobo_FYIProgram_FunctionNet, [351](#)
- GetValve2
 - CRobo_FYIProgram_FunctionNet, [351](#)
- GetValveActive
 - CPPCFunctionNet, [328](#)
 - CWarnerValveControllerDeviceNet, [540](#)
- GetValveActiveEvent
 - CWarnerValveControllerDeviceNet, [553](#)
- GetValveBoardRevision
 - CWarnerValveControllerDeviceNet, [540](#)
- GetValveBoardRevisionEvent
 - CWarnerValveControllerDeviceNet, [553](#)
- GetValveBoardRevisionString
 - CWarnerValveControllerDeviceNet, [540](#)
- GetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, [540](#)
- GetValveDigitalInPortEvent
 - CWarnerValveControllerDeviceNet, [553](#)
- GetValveManualGroup
 - CWarnerValveControllerDeviceNet, [541](#)
- GetValveManualGroupEvent
 - CWarnerValveControllerDeviceNet, [554](#)
- GetValveManualState
 - CWarnerValveControllerDeviceNet, [541](#)
- GetValveManualStateEvent
 - CWarnerValveControllerDeviceNet, [554](#)
- GetValveMode
 - CWarnerValveControllerDeviceNet, [541](#)
- GetValveModeEvent
 - CWarnerValveControllerDeviceNet, [554](#)
- GetValves
 - CMcsBus_FYIExtensionNet, [118](#)
- GetValvesActiveMap
 - CWarnerValveControllerDeviceNet, [542](#)
- GetValvesManualStateMap
 - CWarnerValveControllerDeviceNet, [542](#)
- GetValveTableEntry
 - CWarnerValveControllerDeviceNet, [542](#)
- GetVDD3I
 - CCMOSMea_FunctionNet, [41](#)
- GetVDDI
 - CCMOSMea_FunctionNet, [41](#)
- GetVersion
 - CMcsUsbNet, [235](#)
- GetVersionInt
 - DriverVersionNet, [570](#), [571](#)
- GetVersionString
 - DriverVersionNet, [571](#)
- GetVMMaxNegativeCurrent
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxNegativeCurrentEeprom
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxNegativeVoltage
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxNegativeVoltageEeprom
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxPositiveCurrent
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxPositiveCurrentEeprom
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxPositiveVoltage
 - CMcsBus_VoltageModeNet, [150](#)
- GetVMMaxPositiveVoltageEeprom
 - CMcsBus_VoltageModeNet, [151](#)
- GetVMOutputOnOff
 - CMcsBus_VoltageModeNet, [151](#)
- GetVMVoltage
 - CMcsBus_VoltageModeNet, [151](#)
- GetVoltage
 - COKuvisionStimulatorDeviceNet, [308](#)
- GetVoltage12V
 - CRoboDeviceNet, [373](#)
- GetVoltage12VLimit
 - CRoboDeviceNet, [373](#)
- GetVoltage5V
 - CRoboDeviceNet, [373](#)
- GetVoltage5VLimit
 - CRoboDeviceNet, [373](#)
- GetVoltageAirvalve
 - CRoboDeviceNet, [373](#)
- GetVoltageAirvalveLimit
 - CRoboDeviceNet, [374](#)
- GetVoltageClampControllerParam_D
 - CWarnerUssingFunctionNet, [526](#)
- GetVoltageClampControllerParam_I
 - CWarnerUssingFunctionNet, [526](#)
- GetVoltageClampControllerParam_P
 - CWarnerUssingFunctionNet, [527](#)
- GetVoltageRangeIndex

- CMcsUsbDacqNet, [187](#)
- GetVoltageRangeInMicroVolt
 - CMcsUsbDacqNet, [187](#)
 - CStg200xBasicNet, [429](#)
 - CStimulusFunctionNet, [465](#)
 - CW2100_StimulatorFunctionNet, [510](#)
- GetVoltageRangeInMilliVolt
 - CMcsUsbDacqNet, [188](#)
- GetVoltageResolutionInMicroVolt
 - CStg200xBasicNet, [429](#)
 - CStimulusFunctionNet, [465](#)
 - CW2100_StimulatorFunctionNet, [511](#)
- GetVoltageRs485A
 - CRoboDeviceNet, [374](#)
- GetVoltageRs485ALimit
 - CRoboDeviceNet, [374](#)
- GetVoltageRs485B
 - CRoboDeviceNet, [374](#)
- GetVoltageRs485BLimit
 - CRoboDeviceNet, [374](#)
- GetVoltageValves
 - CRoboDeviceNet, [374](#)
- GetVoltageValvesLimit
 - CRoboDeviceNet, [374](#)
- GetVolti
 - CTcxDeviceNet, [484](#)
- GetWaveform
 - CTEERFunctionNet, [493](#)
- GetWaveLengthInNanometer
 - CMultiwellOptoStimFunctionNet, [300](#)
- GetWorkingFrequency
 - CRFFFunctionNet, [349](#)
- GetWPADebugMode
 - CWClassicFunctionNet, [560](#)
- GetWPAType
 - CWClassicFunctionNet, [560](#)
- GetXGain
 - CRoboDacqNet, [362](#)
- GetXilinxFlashOffset
 - CMcsUsbFactoryNet, [212](#)
- GetXilinxFlashReadCommand
 - CMcsUsbFactoryNet, [212](#)
- GND_SWITCH_BIT
 - CW2100_StimulatorFunctionNet, [513](#)
- GroupID
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [55](#)
- GroupType
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [55](#)
- HasAnalogOut
 - CSCUFunctionNet, [407](#)
- HasGalvanicIsolation
 - CSCUFunctionNet, [407](#)
- HasHSPowerSwitch
 - CSCUFunctionNet, [407](#)
- HasIMU
 - HeadStageIDType, [576](#)
- HasOptoCurrentMeasurement
 - HeadStageIDType, [576](#)
- HasRadioControl
 - CRadioControlledDevicesNet, [342](#)
- HasRef
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- HasRefl
 - CRoboStatorDeviceNet, [389](#)
- HasRefXY
 - CRoboStatorDeviceNet, [389](#)
- HasRefZ
 - CRoboStatorDeviceNet, [389](#)
- HasSoftwareKey
 - CMcsUsbNet, [235](#)
- HeadStageIDType, [575](#)
 - CompareTo, [576](#)
 - ElectricalStimulation, [576](#)
 - Entry, [576](#)
 - Equals, [576](#)
 - HasIMU, [576](#)
 - HasOptoCurrentMeasurement, [576](#)
 - HeadStageIDType, [576](#)
 - HeadstageType, [577](#)
 - HeadstageTypeEnum, [575](#)
 - ID, [577](#)
 - MeasuringOnly, [576](#)
 - NumberOfAnalogChannels, [577](#)
 - NumberOfStimulationChannels, [577](#)
 - OpticalStimulation, [576](#)
 - SN, [577](#)
 - StimulusParameters, [577](#)
 - ToString, [576](#)
 - Type, [577](#)
 - TypeValue, [577](#)
 - Unknown, [576](#)
 - UserDefinedName, [577](#)
 - Valid, [577](#)
 - W16IsW14, [577](#)
- HeadstageIDTypeObject, [578](#)
 - _AdditionalText, [578](#)
 - _IdType, [579](#)
 - AdditionalText, [579](#)
 - Equals, [578](#)
 - GetHashCode, [578](#)
 - HeadstageIDTypeObject, [578](#)
 - IdType, [579](#)
 - ToString, [578](#)
- HeadStageIDTypeState, [579](#)
 - ControlState, [579](#)
 - DataState, [579](#)
 - IdType, [579](#)
 - State, [579](#)
- HeadstageType
 - HeadStageIDType, [577](#)
- HeadstageTypeEnum
 - HeadStageIDType, [575](#)
- HLADacq

- CHLADeviceNet, 101
- HWInfo
 - CMcsUsbDacqNet, 188
- HwVersion
 - CMcsUsbListEntryNet, 220
- ID
 - HeadStageIDType, 577
- IdProduct
 - DeviceIdNet, 565
- IdType
 - HeadstageIDTypeObject, 579
 - HeadstageIDTypeState, 579
- IdVendor
 - DeviceIdNet, 565
- IsAnalogOutEnabled
 - CSCUFunctionNet, 407
- IsAutomaticAnalogOut
 - CSCUFunctionNet, 408
- IsBusy
 - CPPCFunctionNet, 328
- IsChamberAvailable
 - CWarnerUssingFunctionNet, 527
- IsChipPowered
 - CCMOSMea_FunctionNet, 41
- IsConnected
 - CMcsUsbNet, 235
- IsDeviceHighSpeed
 - CMcsUsbNet, 235
- IsDeviceHighSpeedCapable
 - CMcsUsbNet, 236
- IsDeviceTypeOf
 - CMcsUsbListNet, 223
- IsDigitalChannelDedicated
 - CMcsUsbDacqNet::CHWInfo, 103
- IsDigitalOutPortInverted
 - CWarnerValveControllerDeviceNet, 542
- IsDigitalOutPortInvertedEvent
 - CWarnerValveControllerDeviceNet, 554
- IsExceptionsEnabled
 - CMcsUsbNet, 236
- IsGateFloating
 - CCMOSMea_FunctionNet, 41
- IsHeadstageAvailable
 - CSCUFunctionNet, 408
- IsHeadstageAvailableEvent
 - CSCUFunctionNet, 412
- IsHighCurrentMode
 - CWarnerUssingFunctionNet, 527
- IsHSPowered
 - CSCUFunctionNet, 408
- IsInDacqLegacyMode
 - CSCUFunctionNet, 408
- IsInternalCalibrationFinished
 - CTEERFunctionNet, 493
 - CWarnerUssingFunctionNet, 528
- IsPlateTypeValid
 - CMultiwellDeviceNet, 295
- IsPlateTypeValidByHeadstage
 - CMultiwellCallbackFunctionNet, 290
 - CMultiwellDeviceNet, 295
- IsPlateTypeValidByHeadstageEvent
 - CMultiwellCallbackFunctionNet, 291
- IsPulseEnabled
 - CWarnerUssingFunctionNet, 528
- IsPumpMotorOn
 - CRoboFluidDeviceNet, 383
- IsRunning
 - CMeaCleanDeviceNet, 253
 - CMeaCoatDeviceNet, 258
- IsSamplingFinished
 - CTEERFunctionNet, 494
- IsUserTriggerEnabled
 - CLIH3DeviceNet, 112
- IsValveDigitalInInverted
 - CWarnerValveControllerDeviceNet, 543
- IsValveDigitalInInvertedEvent
 - CWarnerValveControllerDeviceNet, 554
- IsValveOpen
 - CWarnerValveControllerDeviceNet, 543
- IsValveOpenEvent
 - CWarnerValveControllerDeviceNet, 554
- IsValveOpenInAnalogMode
 - CWarnerValveControllerDeviceNet, 543
- IsValveOpenInAnalogModeEvent
 - CWarnerValveControllerDeviceNet, 554
- IsValveOpenInDigitalMode
 - CWarnerValveControllerDeviceNet, 544
- IsValveOpenInDigitalModeEvent
 - CWarnerValveControllerDeviceNet, 555
- ListModeSendStart
 - CStg200xBasicNet, 430
- ListModeSendStop
 - CStg200xBasicNet, 430
- ListOfChangedTriggers
 - StgStatusNet, 592
- LoadPressure
 - CPPCFunctionNet, 328
- LoadUserFirmware
 - CMcsUsbFactoryNet, 212
- LoadValveTable
 - CWarnerValveControllerDeviceNet, 544
- LockPlateClamp
 - CMultiwellDeviceNet, 295
- m_Bottom
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, 343
- m_Left
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, 343
- m_pGilsonDevice
 - CGilsonDeviceNet, 98
- m_pMcsBus_MotorControlNet
 - CRoboFluidDeviceNet, 384
- m_pMcsUsb
 - CMcsUsbFunctionNet, 217

- m_pMcsUsbFunction
 - CMcsUsbFunctionNet, [217](#)
- m_pRoboFluidDevice
 - CRoboFluidDeviceNet, [384](#)
- m_Right
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [343](#)
- m_Top
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [343](#)
- Manufacturer
 - CMcsUsbListEntryNet, [220](#)
- MaxBitNumber
 - DigitalSource< digitalsourceenum >, [565](#)
- Mcs, [21](#)
- Mcs::Usb, [21](#)
 - DEVICE_NOT_FOUND, [25](#)
 - enCMosMeaChipType, [25](#)
 - EnSTG200x_STATUS, [25](#)
 - nMos16LV, [25](#)
 - nMos32LV, [25](#)
 - nMos36LN, [25](#)
 - nMos64LN, [25](#)
 - NOT_CONNECTED, [25](#)
 - OK, [25](#)
 - OnChannelData, [25](#)
 - OnDeviceArrivalRemoval, [26](#)
 - OnError, [26](#)
 - OnMcsUsbDeviceState, [26](#)
 - OnMcsUsbDeviceStateCallback, [26](#)
 - OnMwPollStatus, [26](#)
 - OnStg200xDataHandler, [26](#)
 - OnStg200xErrorHandler, [26](#)
 - OnStgPollStatus, [26](#)
 - OnUpdateFirmwareProgress, [27](#)
 - OnUpdateFirmwareStatusChange, [27](#)
 - RoboStatusEventDelegate, [27](#)
 - unknown, [25](#)
- McsBus
 - CPPCDeviceNet, [323](#)
 - CPPS_DeviceNet, [331](#)
 - CRoboDeviceNet, [381](#)
- MCSBUS1
 - FirmwareDestinationNames, [573](#)
- MCSBUS10
 - FirmwareDestinationNames, [573](#)
- MCSBUS11
 - FirmwareDestinationNames, [573](#)
- MCSBUS12
 - FirmwareDestinationNames, [573](#)
- MCSBUS13
 - FirmwareDestinationNames, [573](#)
- MCSBUS2
 - FirmwareDestinationNames, [573](#)
- MCSBUS3
 - FirmwareDestinationNames, [573](#)
- MCSBUS4
 - FirmwareDestinationNames, [574](#)
- MCSBUS5
 - FirmwareDestinationNames, [574](#)
- MCSBUS6
 - FirmwareDestinationNames, [574](#)
- MCSBUS7
 - FirmwareDestinationNames, [574](#)
- MCSBUS8
 - FirmwareDestinationNames, [574](#)
- MCSBUS9
 - FirmwareDestinationNames, [574](#)
- McsBus_MotorControl
 - CPeristalticPumpDeviceNet, [314](#)
 - CPPCDeviceNet, [323](#)
 - CPPS_DeviceNet, [331](#)
 - CRoboDeviceNet, [381](#)
 - CRoboFluidDeviceNet, [384](#)
- McsBus_Sensor
 - CPPCDeviceNet, [323](#)
 - CPPS_DeviceNet, [331](#)
- McsBus_VoltageMode
 - CFluidControlDeviceNet, [78](#)
- McsBus_XY
 - CRoboDeviceNet, [378](#)
- McsBus_ZI
 - CRoboDeviceNet, [378](#)
- McsUsbDeviceStateEvent
 - CMcsUsbDeviceStatePushFunctionNet, [206](#)
 - CMcsUsbDeviceStatePushNet, [207](#)
- MCU1
 - FirmwareDestinationNames, [574](#)
- MeaAudioFunctionNet
 - CMeaDeviceNet, [267](#)
- MeaDigitalDataFunctionNet
 - CMeaDeviceNet, [267](#)
- MeaFeedbackFunctionNet
 - CMeaDeviceNet, [267](#)
- Measure
 - CPathIdentDeviceNet, [311](#)
- MeasureReservoir
 - CPPCFunctionNet, [329](#)
- MeasuringOnly
 - HeadStageIDType, [576](#)
- MeFunctionNet
 - CMeaDeviceNet, [268](#)
- mkfilter
 - mkfilterNet, [580](#)
 - mkfilter_coef_in_one_set
 - mkfilterNet, [580](#)
 - mkfilter_highpass_coeff
 - mkfilterNet, [581](#)
 - mkfilter_highpass_frequency_from_coeff
 - mkfilterNet, [581](#)
 - mkfilter_highpass_frequency_from_k
 - mkfilterNet, [581](#)
 - mkfilter_highpass_k
 - mkfilterNet, [581](#)
 - mkfilter_MCS
 - mkfilterNet, [581](#)

- mkfilter_MCS_k
 - mkfilterNet, [581](#), [582](#)
- mkfilter_normalize_coeffs_int
 - mkfilterNet, [582](#)
- mkfilter_normalize_coeffs_short
 - mkfilterNet, [582](#)
- mkfilter_normalize_scale_coeffs_int
 - mkfilterNet, [582](#)
- mkfilter_scale_coef_in_one_set
 - mkfilterNet, [582](#)
- mkfilterNet, [580](#)
 - mkfilter, [580](#)
 - mkfilter_coef_in_one_set, [580](#)
 - mkfilter_highpass_coeff, [581](#)
 - mkfilter_highpass_frequency_from_coeff, [581](#)
 - mkfilter_highpass_frequency_from_k, [581](#)
 - mkfilter_highpass_k, [581](#)
 - mkfilter_MCS, [581](#)
 - mkfilter_MCS_k, [581](#), [582](#)
 - mkfilter_normalize_coeffs_int, [582](#)
 - mkfilter_normalize_coeffs_short, [582](#)
 - mkfilter_normalize_scale_coeffs_int, [582](#)
 - mkfilter_scale_coef_in_one_set, [582](#)
- MoveAbs
 - CRoboDeviceNet, [374](#)
- MoveAbsI
 - CRoboStatorDeviceNet, [389](#), [390](#)
- MoveAbsXY
 - CRoboStatorDeviceNet, [390](#)
- MoveAbsZ
 - CRoboStatorDeviceNet, [390](#)
- MultibootGetCypressImageld
 - CMcsUsbNet, [236](#)
- MultibootGetImageld
 - CMcsUsbNet, [236](#)
- MultibootGetSelectedImage
 - CMcsUsbNet, [236](#)
- MultibootSelectImage
 - CMcsUsbNet, [236](#)
- MwPollStatusEvent
 - CStg200xDownloadNet, [459](#)
- nMos16LV
 - Mcs::Usb, [25](#)
- nMos32LV
 - Mcs::Usb, [25](#)
- nMos36LN
 - Mcs::Usb, [25](#)
- nMos64LN
 - Mcs::Usb, [25](#)
- NOT_CONNECTED
 - Mcs::Usb, [25](#)
- NumberOfAnalogChannels
 - HeadStageIDType, [577](#)
- NumberOfChannels
 - CDeviceGroupChannelInfoTemplateNet < Dacq-GroupChannelEnumTemplateNet >, [56](#)
- NumberOfStimulationChannels
 - HeadStageIDType, [577](#)
- NumCoefSets
 - CCreateFilterNet, [47](#)
- OK
 - Mcs::Usb, [25](#)
- OnChannelData
 - Mcs::Usb, [25](#)
- OnDeviceArrivalRemoval
 - Mcs::Usb, [26](#)
- OnError
 - Mcs::Usb, [26](#)
- OnGetActiveRunningTableNumber
 - CWarnerValveControllerDeviceNet, [544](#)
- OnGetAnalogThresholdHigh
 - CWarnerValveControllerDeviceNet, [544](#)
- OnGetAnalogThresholdLow
 - CWarnerValveControllerDeviceNet, [544](#)
- OnGetAnalogVoltage
 - CWarnerValveControllerDeviceNet, [544](#)
- OnGetAvailableHeadstages
 - CSCUFunctionNet, [409](#)
- OnGetCurrentNumberOfValves
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetDigitalOutPortValve
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetDisplayMode
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetPlateClampStateByHeadstage
 - CMultiwellCallbackFunctionNet, [290](#)
- OnGetPlateTypeByHeadstage
 - CMultiwellCallbackFunctionNet, [290](#)
- OnGetTableNamebyIndex
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetValveActive
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetValveBoardRevision
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetValveManualGroup
 - CWarnerValveControllerDeviceNet, [545](#)
- OnGetValveManualState
 - CWarnerValveControllerDeviceNet, [546](#)
- OnGetValveMode
 - CWarnerValveControllerDeviceNet, [546](#)
- OnIsDigitalOutPortInverted
 - CWarnerValveControllerDeviceNet, [546](#)
- OnIsHeadstageAvailable
 - CSCUFunctionNet, [409](#)
- OnIsPlateTypeValidByHeadstage
 - CMultiwellCallbackFunctionNet, [290](#)
- OnIsValveDigitalInInverted
 - CWarnerValveControllerDeviceNet, [546](#)
- OnIsValveOpen
 - CWarnerValveControllerDeviceNet, [546](#)
- OnIsValveOpenInAnalogMode
 - CWarnerValveControllerDeviceNet, [546](#)

- OnIsValveOpenInDigitalMode
 - CWarnerValveControllerDeviceNet, [546](#)
- OnMcsUsbDeviceState
 - Mcs::Usb, [26](#)
- OnMcsUsbDeviceStateCallback
 - Mcs::Usb, [26](#)
- OnMwPollStatus
 - Mcs::Usb, [26](#)
- OnStg200xDataHandler
 - Mcs::Usb, [26](#)
- OnStg200xErrorHandler
 - Mcs::Usb, [26](#)
- OnStgPollStatus
 - Mcs::Usb, [26](#)
- OnTableEntryChanged
 - CWarnerValveControllerDeviceNet, [546](#)
- OnUpdateFirmwareProgress
 - Mcs::Usb, [27](#)
- OnUpdateFirmwareStatusChange
 - Mcs::Usb, [27](#)
- OpenPipe
 - CGenericDevelopDeviceNet, [91](#)
- OpenPlateClamp
 - CMultiwellDeviceNet, [295](#)
- operator=
 - DeviceIdNet, [564](#)
- OpticalStimulation
 - HeadStageIDType, [576](#)
- Order
 - CCreateFilterNet, [47](#)
 - CFilterPropertyNet, [71](#)
- PatternListStart
 - COctoPotDeviceNet, [304](#)
- PIC
 - FirmwareDestinationNames, [574](#)
- PIC2
 - FirmwareDestinationNames, [574](#)
- PIC3
 - FirmwareDestinationNames, [574](#)
- PIC4
 - FirmwareDestinationNames, [574](#)
- PollStatusEvent
 - CStimulusFunctionNet, [471](#)
 - CW2100_StimulatorFunctionNet, [513](#)
- PowerChip
 - CCMOSMea_FunctionNet, [41](#)
- PowerHS
 - CSCUFunctionNet, [409](#)
- PPCFunction
 - CPPCDeviceNet, [323](#)
- PPS_Function
 - CPPS_DeviceNet, [331](#)
- PrepareAndAppendData
 - CStg200xDownloadNet, [456](#)
 - CStimulusFunctionNet, [466](#)
- PrepareAndSendData
 - CStg200xDownloadNet, [457](#)
 - CStimulusFunctionNet, [467](#)
- PrepareChannelData
 - CDigOutStimulatorFunctionNet, [59](#)
- PrepareData
 - CStimulusFunctionNet, [467](#)
 - CW2100_StimulatorFunctionNet, [511](#)
- PrepareDataSync
 - CW2100_StimulatorFunctionNet, [511](#)
- Product
 - CMcsUsbListEntryNet, [221](#)
- Program
 - CProgramPressureCurveNet, [337](#)
- PulseGenerator
 - CW2100_FunctionNet, [506](#)
- PumpOff
 - CRoboFluidDeviceNet, [383](#)
- PumpOn
 - CRoboFluidDeviceNet, [383](#)
- QueryTriggerstatus
 - CStg200xDownloadNet, [458](#)
- RampStart
 - COctoPotDeviceNet, [304](#)
- Read
 - CExternDTesterDeviceNet, [62](#)
- Read2
 - CExternDTesterDeviceNet, [63](#)
- ReadBlockFromFlash
 - CMcsUsbFactoryNet, [212](#)
- ReadBlockFromIFBGlobalEeprom
 - CMcsUsbFactoryNet, [213](#)
- ReadBlockFromNVMEM
 - CMcsUsbFactoryNet, [213](#)
- ReadClipping
 - CLIH3DeviceNet, [112](#)
- ReadEepromRegisterPreconfig
 - CMcsUsbNet, [237](#)
- ReadPipe
 - CGenericDevelopDeviceNet, [91](#)
- ReadRegister
 - CMcsUsbNet, [237](#)
- ReadRegister32
 - CMcsUsbNet, [237](#)
- ReadRegisterTimeSlot
 - CMcsUsbNet, [237](#)
- Receive
 - CSerialPortNet, [413](#)
- ReceiveString
 - CSerialPortNet, [413](#)
- RemoveSoftwareKey
 - CMcsUsbNet, [238](#)
- RescanHeadstage
 - CMcsUsbNet, [238](#)
- ResetAdcOffset
 - COctoPotDeviceNet, [304](#)
- ResetChannelmap
 - CWClassicFunctionNet, [561](#)
- ResetDacOffset
 - COctoPotDeviceNet, [304](#)

- ResetHighpassFilter
 - CFilterConfigurationNet, [67](#)
- ResetPipe
 - CGenericDevelopDeviceNet, [91](#)
- ResetStatus
 - CStg200xDownloadBasicNet, [448](#)
- RFFunction
 - CPositionIIDeviceNet, [319](#)
- RoboDacq
 - CHiClampDeviceNet, [99](#)
- RoboDevice
 - CSafeISDeviceNet, [396](#)
- RoboError_AnotherMaster
 - CRoboDeviceNet, [378](#)
- RoboError_Base
 - CRoboDeviceNet, [378](#)
- RoboError_CannotEscapeEndSwitch
 - CRoboDeviceNet, [378](#)
- RoboError_CommandAlreadyInProgress
 - CRoboDeviceNet, [378](#)
- RoboError_CommandNotPossible
 - CRoboDeviceNet, [378](#)
- RoboError_CommunicationTimeout
 - CRoboDeviceNet, [378](#)
- RoboError_DacqNotReady
 - CRoboDeviceNet, [379](#)
- RoboError_DLLMovementTimeout
 - CRoboDeviceNet, [379](#)
- RoboError_FindReferenceMethod
 - CRoboDeviceNet, [379](#)
- RoboError_GilsonCommandPending
 - CRoboDeviceNet, [379](#)
- RoboError_GilsonTimeout
 - CRoboDeviceNet, [379](#)
- RoboError_GilsonWrondID
 - CRoboDeviceNet, [379](#)
- RoboError_McsBus_UnknownCommand
 - CRoboDeviceNet, [379](#)
- RoboError_NoEndSwitch
 - CRoboDeviceNet, [379](#)
- RoboError_NoMoreData
 - CRoboDeviceNet, [379](#)
- RoboError_NoReference
 - CRoboDeviceNet, [379](#)
- RoboError_NoSpeedOrAcceleration
 - CRoboDeviceNet, [380](#)
- RoboError_OverPressure
 - CRoboDeviceNet, [380](#)
- RoboError_ParameterNotAllowed
 - CRoboDeviceNet, [380](#)
- RoboError_PeristalticTimeout
 - CRoboDeviceNet, [380](#)
- RoboError_Phase0OutOfRange
 - CRoboDeviceNet, [380](#)
- RoboError_PollLoopCanceled
 - CRoboDeviceNet, [380](#)
- RoboError_PollLoopCanceledAndStopMovement
 - CRoboDeviceNet, [380](#)
- RoboError_Pressure
 - CRoboDeviceNet, [380](#)
- RoboError_RangeExceeded
 - CRoboDeviceNet, [380](#)
- RoboError_StateChangeNotPossible
 - CRoboDeviceNet, [380](#)
- RoboError_Timeout
 - CRoboDeviceNet, [381](#)
- RoboError_UnknownCommand
 - CRoboDeviceNet, [381](#)
- RoboMainLowLevelCommand
 - CRoboDeviceNet, [381](#)
- RoboMainStatorLowLevelCommand
 - CRoboStatorDeviceNet, [393](#)
- RoboStatusEvent
 - CRoboDeviceNet, [381](#)
- RoboStatusEventDelegate
 - Mcs::Usb, [27](#)
- RunTable
 - CRoboDacqNet, [362](#)
- SampleRate
 - CCreateFilterNet, [47](#)
- Samplerate
 - CMcsUsbDacqNet, [204](#)
- Scale
 - CCreateFilterNet, [47](#)
- ScanForHeadstages
 - CWClassicFunctionNet, [561](#)
- SelectHeadstage
 - CW2100_FunctionNet, [505](#)
- SelectTimeSlot
 - CW2100_StimulatorFunctionNet, [511](#)
- Send
 - CSerialPortNet, [413](#)
- SendBuffered
 - CGilsonDeviceNet, [98](#)
- SendChannelData
 - CDigOutStimulatorFunctionNet, [59](#)
 - CStg200xDownloadBasicNet, [450](#)
- SendCommand
 - CLIH3DeviceNet, [112](#)
- SendImmediate
 - CGilsonDeviceNet, [98](#)
- SendImmediateGetResponse
 - CGilsonDeviceNet, [98](#)
- SendMultiplexedData
 - CStimulusFunctionNet, [468](#)
- SendPreparedData
 - CStimulusFunctionNet, [468](#)
 - CW2100_StimulatorFunctionNet, [511](#)
- SendSegmentDefine
 - CStg200xDownloadNet, [458](#)
- SendSegmentSelect
 - CStg200xDownloadNet, [458](#)
- SendSegmentStart
 - CStg200xDownloadNet, [459](#)
- SendStart
 - CStg200xBasicNet, [430](#)

- CStimulusFunctionNet, [468](#)
- CW2100_StimulatorFunctionNet, [512](#)
- SendStartDacq
 - CMcsUsbDacqNet, [188](#)
- SendStartStgAndDacq
 - CMcsUsbDacqNet, [188](#)
- SendStop
 - CStg200xBasicNet, [430](#)
 - CStimulusFunctionNet, [468](#)
 - CW2100_StimulatorFunctionNet, [512](#)
- SendStopDacq
 - CMcsUsbDacqNet, [189](#)
- SendStopStgAndDacq
 - CMcsUsbDacqNet, [189](#)
- SendStopStgAndDacqWithOptions
 - CMcsUsbDacqNet, [189](#)
- SendSyncData
 - CStg200xDownloadBasicNet, [450](#)
- Sensor
 - CFYIDeviceNet, [79](#)
 - CMeasureTableDeviceNet, [276](#)
 - CPatchServerDeviceNet, [310](#)
- SerialNumber
 - CMcsUsbListEntryNet, [221](#)
 - CMcsUsbNet, [245](#)
- SerialPort
 - CHLADeviceNet, [101](#)
- Set4ADCCatchampAverageShift
 - CMcsBus_SensorNet, [144](#)
- Set4ADCMode
 - CMcsBus_SensorNet, [144](#)
- Set4DAC
 - CMcsBus_SensorNet, [144](#)
- Set_Values
 - CNF_GenDeviceNet, [302](#)
 - CPathIdentDeviceNet, [311](#)
- SetAbsMaxCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [300](#)
- SetAccelerationI
 - CRoboStatorDeviceNet, [390](#)
- SetAccelerationNativeI
 - CRoboStatorDeviceNet, [390](#)
- SetAccelerationNativeXY
 - CRoboStatorDeviceNet, [390](#)
- SetAccelerationNativeZ
 - CRoboStatorDeviceNet, [390](#)
- SetAccelerationXY
 - CRoboStatorDeviceNet, [391](#)
- SetAccelerationZ
 - CRoboStatorDeviceNet, [391](#)
- SetAccelGyroDesiredRate
 - CW2100_FunctionNet, [505](#)
- SetAccelGyroEnabled
 - CW2100_FunctionNet, [505](#)
- SetAccelRange
 - CW2100_FunctionNet, [505](#)
- SetActiveRunningTableNumber
 - CWarnerValveControllerDeviceNet, [546](#)
- SetADC
 - CWarnerValveControllerDeviceTesterFunctionNet, [557](#)
- SetAdcChannels
 - CSafeISDeviceNet, [394](#)
- SetADCInputOffset
 - CCMOSMea_FunctionNet, [41](#)
- SetADCOffset
 - CLIH3DeviceNet, [113](#)
- SetAdcOffset
 - COctoPotDeviceNet, [305](#)
- SetAdcSamplePos
 - CSafeISDeviceNet, [394](#)
- SetAFormat
 - CFilterCoefficientsNet, [65](#)
- SetAirpressureLimit
 - CRoboDeviceNet, [375](#)
- SetAirValve
 - CRoboDeviceNet, [375](#)
- SetAllDigout
 - CRoboDacqNet, [363](#)
- SetAmplificationSwitch
 - COctoPotDeviceNet, [305](#)
- SetAmplitude
 - CChannelTestDeviceNet, [31](#)
- SetAmplitude_nA
 - CTEERFunctionNet, [494](#)
- SetAnalogOutADCRange
 - CSCUFunctionNet, [409](#)
- SetAnalogOutChannel
 - CW2100_FunctionNet, [505](#)
- SetAnalogOutChannels
 - CSCUFunctionNet, [409](#)
- SetAnalogOutDACRange
 - CSCUFunctionNet, [411](#)
- SetAnalogOutFilter
 - CW2100_FunctionNet, [505](#)
- SetAnalogThresholdHigh
 - CWarnerValveControllerDeviceNet, [547](#)
- SetAnalogThresholdLow
 - CWarnerValveControllerDeviceNet, [547](#)
- SetAnalogVoltageRange
 - CPPCFunctionNet, [329](#)
- SetAnalogVoltages
 - CPPS_FunctionNet, [334](#)
- SetAttenuation
 - CChannelTestDeviceNet, [31](#)
- SetAudioChannels
 - CMeaAudioFunctionNet, [249](#), [250](#)
 - CW2100_FunctionNet, [505](#)
- SetAutocalibrationDisabled
 - CStg200xBasicNet, [430](#)
- SetAxisConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetAxisLED
 - CRoboocyte2DeviceNet, [387](#)
- SetAxisParametersEeprom

- CMcsBus_AxisParametersNet, [115](#), [116](#)
- SetBandwidthByIndex
 - CIantanMea_FunctionNet, [105](#)
- SetBaseFrequency
 - CRFFunctionNet, [349](#)
- SetBaseSamplerate
 - CCMOSMeaDeviceNet, [45](#)
- SetBath
 - CCMOSMea_FunctionNet, [41](#)
- SetBathclamp
 - COctoPotDeviceNet, [305](#)
- SetBathMode
 - CCMOSMea_FunctionNet, [41](#)
- SetBFormat
 - CFilterCoefficientsNet, [65](#)
- SetBlankingEnable
 - CStg200xBasicNet, [431](#), [432](#)
- SetBoostAlwaysOnMode
 - CW2100_StimulatorFunctionNet, [512](#)
- SetBuffer
 - CGenericDevelopDeviceNet, [92](#)
- SetBufferIndex
 - CTEERFunctionNet, [494](#)
- SetBusAddress
 - CMcsBusNet, [155](#)
- SetBusAddressEeprom
 - CMcsBusNet, [155](#)
- SetByteBuffer
 - CGenericDevelopDeviceNet, [92](#)
- SetCalibration
 - CTcxDeviceNet, [484](#)
- SetCardinalDacqSamplerate
 - CInterfaceboardFunctionNet, [107](#)
- SetCardinalStgOutputrate
 - CInterfaceboardFunctionNet, [107](#)
- SetChannel
 - CSw2to64DeviceNet, [473](#)
- SetChannelmap
 - CWClassicFunctionNet, [561](#)
- SetChannels
 - CSw2to64DeviceNet, [473](#)
- SetChannelSwitch
 - COctoPotDeviceNet, [305](#)
- SetChargingMode
 - CMultiBatteryChargerDeviceNet, [286](#)
- SetChargingPCoefficient
 - CMultiBatteryChargerDeviceNet, [286](#)
- SetCheckVoltage
 - CokuvisionStimulatorDeviceNet, [309](#)
- SetClampMode
 - CTEERFunctionNet, [494](#)
 - CWarnerUssingFunctionNet, [528](#)
- SetColorRgb
 - CMultiwellOptoStimFunctionNet, [300](#)
- SetColorStr
 - CMultiwellOptoStimFunctionNet, [300](#)
- SetCommand
 - CMcsBusNet, [156](#)
- CPedoterDeviceNet, [313](#)
- CRoboDacqNet, [363](#)
- SetConfiguration
 - CMcsUsbNet, [238](#)
- SetConfigurationBit
 - CRoboDacqNet, [363](#)
- SetConfigurationBitAxc
 - CRoboDacqNet, [363](#)
- SetConfigurationBitBlu_Led
 - CRoboDacqNet, [363](#)
- SetConfigurationBitBlu_LedToggleFast
 - CRoboDacqNet, [363](#)
- SetConfigurationBitBlu_LedToggleSlow
 - CRoboDacqNet, [363](#)
- SetConfigurationBitCC_Gen
 - CRoboDacqNet, [363](#)
- SetConfigurationBitCV_Gen
 - CRoboDacqNet, [363](#)
- SetConfigurationBitRC_Gen
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRed_Led
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRed_LedSaturation
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRed_LedToggleFast
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRed_LedToggleSlow
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRelais
 - CRoboDacqNet, [364](#)
- SetConfigurationBitRV_Gen
 - CRoboDacqNet, [364](#)
- SetConfigurationBitStream
 - CRoboDacqNet, [364](#)
- SetConfigurationBitSupply
 - CRoboDacqNet, [364](#)
- SetControllerParams
 - CTEERFunctionNet, [495](#)
- SetCrossTalkOffset
 - CRoboDacqNet, [364](#)
- SetCrossTalkOptimum
 - CRoboDacqNet, [365](#)
- SetCurrentAirvalveLimit
 - CRoboDeviceNet, [375](#)
- SetCurrentAndAir
 - CRoboDeviceNet, [375](#)
- SetCurrentAndAirXY
 - CRoboStatorDeviceNet, [391](#)
- SetCurrentEditTableNumber
 - CWarnerValveControllerDeviceNet, [547](#)
- SetCurrentEnable
 - CTEERFunctionNet, [495](#)
- SetCurrentFactor
 - CokuvisionStimulatorDeviceNet, [309](#)
- SetCurrentMode
 - CStg200xBasicNet, [432](#)
- SetCycles
 - CMeaCleanDeviceNet, [253](#)

- CMeaCoatDeviceNet, 258
- SetD
 - CTcxDeviceNet, 485
- SetDacAmplificationFactor
 - CStg200xBasicNet, 432
- SetDacAutoControl
 - COctoPotDeviceNet, 305
- SetDacIdleValue
 - CLIH3DeviceNet, 113
- SetDacMode
 - CSafeISDeviceNet, 394
- SetDACOffset
 - CokuvisionStimulatorDeviceNet, 309
- SetDacOffset
 - CDacCalibrationFunctionNet, 49
 - COctoPotDeviceNet, 305
- SetDacPeriode
 - CSafeISDeviceNet, 395
- SetDacPulseform
 - CSafeISDeviceNet, 395
- SetDacqLegacyMode
 - CSCUFunctionNet, 411
- SetDacRange
 - CW2100_FunctionNet, 506
- SetDACs
 - CMcsBus_SensorNet, 144
- SetDacUseIdleValue
 - CLIH3DeviceNet, 113
- SetDacValue
 - COctoPotDeviceNet, 305
- SetDataMode
 - CMcsUsbDacqNet, 190
- SetDefault
 - CWarnerValveControllerDeviceNet, 547
- SetDestinationSerialNumber
 - CMcsUsbFactoryNet, 213
- SetDetectionThreshold
 - CMcsBus_SensorNet, 144
- SetDevice
 - CTcxDeviceNet, 485
- SetDeviceList
 - CPositionImpDeviceNet, 322
- SetDeviceType
 - CTcxDeviceNet, 485
- SetDevname
 - CTcxDeviceNet, 485
- SetDiagnosticMode
 - CIntanMea_FunctionNet, 105
- SetDigitalData
 - CMeaDigitalDataFunctionNet, 269
- SetDigitalOut
 - CMeaDeviceNet, 264
- SetDigitalOutPortInvert
 - CWarnerValveControllerDeviceNet, 547
- SetDigitalOutPortValve
 - CWarnerValveControllerDeviceNet, 548
- SetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, 548
- SetDigitalSource
 - CMcsUsbDacqNet, 190–192
- SetDigitalStimulatorTrigger
 - CW2100_StimulatorFunctionNet, 512
- SetDigitalStimulatorTriggerSlope
 - CW2100_StimulatorFunctionNet, 512
- SetDigout
 - CFluidControlDeviceNet, 76
 - CRoboDacqNet, 365
- SetDigoutMode
 - CStg200xBasicNet, 433
- SetDigOutState
 - CLIH3DeviceNet, 113
- SetDigoutValue
 - CStg200xBasicNet, 433
- SetDIO
 - CMcsBus_FYIExtensionNet, 118
- SetDischargeCurrentSetPoint
 - CMultiBatteryChargerDeviceNet, 287
- SetDisplayMode
 - CWarnerValveControllerDeviceNet, 548
- SetDisplayText
 - CRoboDacqNet, 365
- SetDownsampleFactor
 - CRoboDacqNet, 365
- SetDSPHighPassByIndex
 - CIntanMea_FunctionNet, 105
- SetDuration
 - CMeaCoatDeviceNet, 258
- SetEEPromPage
 - CLIH3DeviceNet, 114
- SetElectrodeDacMux
 - CStg200xBasicNet, 433–435
- SetElectrodeEnable
 - CStg200xBasicNet, 435–437
- SetElectrodeMode
 - CStg200xBasicNet, 438, 439
- SetEnableAmplifierProtectionSwitch
 - CStg200xBasicNet, 439, 440
- SetEnableHeaterLimit
 - CTcxDeviceNet, 485
- SetEnablePulse
 - CWarnerUssingFunctionNet, 529
- SetEnableThermocouple
 - CTcxDeviceNet, 485
- SetExternalElectrodeEnable
 - CStg200xBasicNet, 441
- SetExternalLED
 - CTEERFunctionNet, 495
- SetFAAmplification
 - CStg200xBasicNet, 442
- SetFilter
 - CRoboDacqNet, 365
- SetFilterCoeffs
 - CRoboDacqNet, 365
- SetFilterParameter
 - CFilterConfigurationNet, 68
 - CFilterConfigurationRegisterNet, 69

- SetFilterParameterPermanent
 - CFilterConfigurationNet, [68](#)
 - CFilterConfigurationRegisterNet, [69](#), [70](#)
- SetFilterParametersHeadstage
 - CWClassicFunctionNet, [561](#)
- SetFinalDischargeVoltage
 - CMultiBatteryChargerDeviceNet, [287](#)
- SetFrequency
 - CChannelTestDeviceNet, [31](#)
 - CRadioControlledDevicesNet, [342](#)
- SetGain
 - CPgaDeviceNet, [316](#)
- SetGate
 - CCMOSMea_FunctionNet, [41](#)
- SetGateFloating
 - CCMOSMea_FunctionNet, [42](#)
- SetGateToVOP
 - CCMOSMea_FunctionNet, [42](#)
- SetGlobalRepeat
 - CDigOutStimulatorFunctionNet, [60](#)
- SetGyroRange
 - CW2100_FunctionNet, [506](#)
- SetHasChecksum
 - CWClassicFunctionNet, [561](#)
- SetHeadstage
 - CStg200xBasicNet, [442](#)
- SetHeadstageOnOff
 - CW2100_FunctionNet, [506](#)
 - CWClassicFunctionNet, [561](#)
- SetHeadstageSamplingActive
 - CW2100_FunctionNet, [506](#)
- SetHeadstageToSleep
 - CW2100_FunctionNet, [506](#)
- SetHeaterLimit
 - CTcxDeviceNet, [485](#)
- SetHighCurrentMode
 - CWarnerUssingFunctionNet, [529](#)
- SetHighpassFilterEnable
 - CFilterConfigurationNet, [68](#)
- SetHWConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetHWRevision
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetHWRevisionEeprom
 - CMcsBusNet, [156](#)
- SetHWSelectedChannels
 - CWClassicFunctionNet, [561](#)
- SetI
 - CTcxDeviceNet, [485](#)
- SetIClamp
 - CRoboDacqNet, [365](#)
- SetICoeff
 - CRobo_FYITemp_FunctionNet, [353](#)
- SetICOffset
 - CRoboDacqNet, [365](#)
- SetIdleModeOffset
 - CWarnerUssingFunctionNet, [529](#)
- SetIGain
 - CRoboDacqNet, [365](#)
- SetImpedanceTestFrequency
 - CMealImpedanceDeviceNet, [275](#)
- SetImpId
 - CPositionImpDeviceNet, [322](#)
- SetImplantCurrentSetpoint
 - CPositionIIDeviceNet, [319](#)
- SetInMovement
 - CRoboDeviceNet, [375](#)
- SetIntanRegister
 - CIntanMea_FunctionNet, [105](#)
- SetIntBuffer
 - CGenericDevelopDeviceNet, [92](#)
- SetIO
 - CWarnerValveControllerDeviceTesterFunctionNet, [557](#)
- SetIODirection
 - CWarnerValveControllerDeviceTesterFunctionNet, [557](#)
- SetLatency
 - CMcsBus_SensorNet, [144](#)
- SetLayoutConfiguration
 - CMEA2100x256FunctionNet, [246](#)
- SetLED
 - CRetinaLedDeviceNet, [345](#)
- SetLEDSwitch
 - CMcsBus_ExtensionNet, [117](#)
- SetLength
 - CRobo_FYIProgram_FunctionNet, [351](#)
- SetLiquidResistance
 - CTEERFunctionNet, [495](#)
 - CWarnerUssingFunctionNet, [529](#)
- SetListmodeIndexRange
 - CStg200xBasicNet, [442](#)
- SetListmodeTriggerSource
 - CStg200xBasicNet, [442](#)
- SetLowCurrentMode
 - CWarnerUssingFunctionNet, [530](#)
- SetLumi
 - CRetinaLedDeviceNet, [345](#)
- SetMaxCurrent
 - CMeaCoatDeviceNet, [259](#)
- SetMaxDurationHighCurrentInMicroSec
 - CMultiwellOptoStimFunctionNet, [301](#)
- SetMaxDutyCycleHighCurrent
 - CMultiwellOptoStimFunctionNet, [301](#)
- SetMaxHeaterPowerMultiwell
 - CTcxDeviceNet, [486](#)
- SetMaxNoPressure
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMaxNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMaxP
 - CTcxDeviceNet, [486](#)

- SetMaxPower
 - CokuvisionStimulatorDeviceNet, [309](#)
 - CRobo_FYITemp_FunctionNet, [353](#)
- SetMaxPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMaxVoltage
 - CMeaCleanDeviceNet, [253](#)
 - CokuvisionStimulatorDeviceNet, [309](#)
- SetMCAcceleration
 - CMcsBus_MotorControlNet, [128](#)
- SetMCAccelerationEeprom
 - CMcsBus_MotorControlNet, [129](#)
- SetMCAccelerationShortCommand
 - CMcsBus_MotorControlNet, [129](#)
- SetMCAXisRevisionEeprom
 - CMcsBus_MotorControlNet, [129](#)
- SetMCBreakCurrent
 - CMcsBus_MotorControlNet, [129](#)
- SetMCBreakCurrentEeprom
 - CMcsBus_MotorControlNet, [129](#)
- SetMCConfig
 - CMcsBus_MotorControlNet, [129](#)
- SetMCConfigEeprom
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrent
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrentEeprom
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrentMode
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrentModeEeprom
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrentModeShortCommand
 - CMcsBus_MotorControlNet, [130](#)
- SetMCCurrentPosition
 - CMcsBus_MotorControlNet, [131](#)
- SetMCCurrentShortCommand
 - CMcsBus_MotorControlNet, [131](#)
- SetMCMaxAcceleration
 - CMcsBus_MotorControlNet, [131](#)
- SetMCMaxAccelerationEeprom
 - CMcsBus_MotorControlNet, [131](#)
- SetMCMaxCurrent
 - CMcsBus_MotorControlNet, [131](#)
- SetMCMaxCurrentEeprom
 - CMcsBus_MotorControlNet, [131](#)
- SetMCMaxSpeed
 - CMcsBus_MotorControlNet, [132](#)
- SetMCMaxSpeedEeprom
 - CMcsBus_MotorControlNet, [132](#)
- SetMCMaxTravel
 - CMcsBus_MotorControlNet, [132](#)
- SetMCMaxTravelEeprom
 - CMcsBus_MotorControlNet, [132](#)
- SetMCMaxTravelShortCommand
 - CMcsBus_MotorControlNet, [132](#)
- SetMCNewPosition
 - CMcsBus_MotorControlNet, [132](#)
- SetMCOutputOnOff
 - CMcsBus_MotorControlNet, [133](#)
- SetMCReference
 - CMcsBus_MotorControlNet, [133](#)
- SetMCReferenceCurrent
 - CMcsBus_MotorControlNet, [133](#)
- SetMCReferenceCurrentEeprom
 - CMcsBus_MotorControlNet, [133](#)
- SetMCRegulatorGain
 - CMcsBus_MotorControlNet, [133](#)
- SetMCRegulatorGainEeprom
 - CMcsBus_MotorControlNet, [133](#)
- SetMCRotation
 - CMcsBus_MotorControlNet, [134](#)
- SetMCScalingFactor
 - CMcsBus_MotorControlNet, [134](#)
- SetMCScalingFactorEeprom
 - CMcsBus_MotorControlNet, [134](#)
- SetMCSpeed
 - CMcsBus_MotorControlNet, [134](#)
- SetMCSpeedEeprom
 - CMcsBus_MotorControlNet, [134](#)
- SetMCSpeedShortCommand
 - CMcsBus_MotorControlNet, [134](#)
- SetMCSpeedUnitEeprom
 - CMcsBus_MotorControlNet, [135](#)
- SetMCStandbyCurrent
 - CMcsBus_MotorControlNet, [135](#)
- SetMCStandbyCurrentEeprom
 - CMcsBus_MotorControlNet, [135](#)
- SetMCStandbyTime
 - CMcsBus_MotorControlNet, [135](#)
- SetMCStandbyTimeEeprom
 - CMcsBus_MotorControlNet, [135](#)
- SetMeasurementMode
 - CStg200xBasicNet, [442](#)
- SetMinimalThreshold
 - CMcsBus_SensorNet, [144](#)
- SetMinNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMinPressure
 - CRoboDeviceNet, [375](#)
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMinPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [587](#)
- SetMinVoltage
 - CMeaCleanDeviceNet, [254](#)
- SetModeSelect
 - CPulseGeneratorFunctionNet, [340](#)
- SetMovePump
 - CMcsBus_SensorNet, [145](#)
- SetMultiHeadstageMode
 - CW2100_FunctionNet, [506](#)
- SetNanoVoltsPerKelvin

- CMcsBus_TempSensorNet, 148
- SetNeurochipMemoryData
 - CCMOSMea_FunctionNet, 42
- SetNumberOfAnalogChannels
 - CMeaDeviceNet, 264
- SetNumberOfChannels
 - CMeaDeviceNet, 265, 266
 - COctoPotDeviceNet, 305
- SetOffsetCurrent
 - CMeaCoatDeviceNet, 259
- SetOnOff
 - CTcxDeviceNet, 486
- SetOutputMap
 - CStg200xDownloadNet, 459
- SetOutputRate
 - COctoPotDeviceNet, 305
 - CStg200xBasicNet, 443
- SetP
 - CTcxDeviceNet, 486
- SetParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetPattern
 - CMeaSwitchDeviceNet, 278
- SetPatternBool
 - CMeaSwitchDeviceNet, 278
- SetPatternListEntry
 - COctoPotDeviceNet, 306
- SetPauseDuration
 - CMeaCoatDeviceNet, 259
- SetPCoeff
 - CRobo_FYITemp_FunctionNet, 353
- SetPeriod
 - CPulseGeneratorFunctionNet, 340
- SetPeriod_us
 - CTEERFunctionNet, 496
- SetPermanentCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, 301
- SetPersistency
 - CRetinaLedDeviceNet, 345
- SetPGain
 - CRoboDacqNet, 366
- SetPidParameter
 - COctoPotDeviceNet, 306
- SetPiezoState
 - CMcsBus_SensorNet, 145
- SetPlateMux
 - CMultiwellDeviceNet, 295
- SetPlateMuxByHeadstage
 - CMultiwellDeviceNet, 296
- SetPlateType
 - CMultiwellDeviceNet, 296
- SetPlateTypeByHeadstage
 - CMultiwellDeviceNet, 296
- SetPoti
 - CMcsUsbDacqNet, 192
- SetPressureOffset
 - CMcsBus_SensorNet, 145
- CPPCFunctionNet, 329
- SetPressureRange
 - CPPCFunctionNet, 329
- SetPulse
 - CWarnerUssingFunctionNet, 530
- SetPulseform
 - COkuvisionStimulatorDeviceNet, 309
- SetPulseLength
 - CPulseGeneratorFunctionNet, 341
- SetPumpCouple
 - CPPS_FunctionNet, 334
- SetPumpEnableSpeedRatio
 - CPPS_FunctionNet, 334
- SetPumpFastOnOff
 - CPPS_FunctionNet, 334
- SetPumpFastSpeed
 - CPPS_FunctionNet, 335
- SetPumpFunctionSpeeds
 - CPPS_FunctionNet, 335
- SetPumpManualOnOff
 - CPPS_FunctionNet, 335
- SetPumpMaxSpeed
 - CPPS_FunctionNet, 335
- SetPumpModeType
 - CPPCFunctionNet, 329
 - CPPS_FunctionNet, 335
- SetPumpSpeed
 - CRoboFluidDeviceNet, 383
- SetPumpSpeedRatio
 - CPPS_FunctionNet, 335
- SetPumpSpeedUnit
 - CPPCFunctionNet, 330
 - CPPS_FunctionNet, 335
- SetPWM
 - CFluidControlDeviceNet, 76
- SetRampParameter
 - COctoPotDeviceNet, 306
- SetRatedCapacity
 - CMultiBatteryChargerDeviceNet, 287
- SetRatedCapacityVolatile
 - CMultiBatteryChargerDeviceNet, 287
- SetRecordingNumber
 - CRoboDacqNet, 366
- SetReferenceElectrodeMode
 - CSCUFunctionNet, 411
- SetReferenceElectrodeSwitchState
 - CSCUFunctionNet, 411
- SetRegionOfInterests
 - CCMOSMeaDeviceNet, 45
- SetRegulationTimeouts
 - CMcsBus_SensorNet, 145
- SetRegulatorFactor
 - CMcsBus_SensorNet, 145
- SetRegulatorOnOff
 - CMcsBus_SensorNet, 145
 - CRobo_FYITemp_FunctionNet, 353
- SetRepeat
 - CRetinaLedDeviceNet, 345

- SetRepeats
 - CProgramPressureCurveNet, 338
- SetResetFilter
 - CWClassicFunctionNet, 561
- SetRFFrequency
 - CPositionImpDeviceNet, 322
- SetRFFrequencyHeadstage
 - CWClassicFunctionNet, 561
- SetRFFrequencyReceiver
 - CWClassicFunctionNet, 562
- SetRFFrequencyReceiverEeprom
 - CWClassicFunctionNet, 562
- SetRFLostBehaviour
 - CWClassicFunctionNet, 562
- SetRFPower
 - CWClassicFunctionNet, 562
- SetRotatePump
 - CMcsBus_SensorNet, 145
- SetRTC
 - CokuvisionStimulatorDeviceNet, 309
- SetSampleInterval
 - CLIH3DeviceNet, 114
- SetSamplePeriode
 - CMcsBus_SensorNet, 146
- SetSamplerate
 - CMcsUsbDacqNet, 192
- SetScreen
 - CRoboDacqNet, 366
- SetSearchReferenceFastAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetSearchReferenceFastSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetSearchReferenceFineAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetSearchReferenceFineSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetSearchReferenceMethod
 - CRoboDeviceNet::RoboMainLowLevelCommands, 588
- SetSearchReferenceMoveOut
 - CRoboDeviceNet::RoboMainLowLevelCommands, 589
- SetSearchReferenceOffsetPos
 - CRoboDeviceNet::RoboMainLowLevelCommands, 589
- SetSelectedChannels
 - CMcsUsbDacqNet, 192–194
 - CW2100_FunctionNet, 506
- SetSelectedChannelsQueue
 - CMcsUsbDacqNet, 195–197
- SetSelectedData
 - CMcsUsbDacqNet, 197–199
- SetSelectedHeadstage
 - CWClassicFunctionNet, 562
- SetSensorType
 - CTcxDeviceNet, 486
- SetSerialNumberHeadstage
 - CWClassicFunctionNet, 562
- SetSetpoint
 - CTcxDeviceNet, 486
- SetShortBuffer
 - CGenericDevelopDeviceNet, 93
- SetSimulation
 - CRoboDacqNet, 366
- SetSineParameter
 - COctoPotDeviceNet, 306
- SetSingleHeater
 - CMcsBus_FYIExtensionNet, 118
- SetSingleValve
 - CFluidControlDeviceNet, 77
 - CRoboFluidDeviceNet, 384
- SetSlope
 - CMeaCleanDeviceNet, 254
 - CMeaCoatDeviceNet, 259
- SetSoftwareKey
 - CMcsUsbNet, 238
- SetSollPressure
 - CMcsBus_SensorNet, 146
- SetSollTemp
 - CRobo_FYITemp_FunctionNet, 354
- SetSourceBulk
 - CCMOSMea_FunctionNet, 42
- SetSourceDrain
 - CCMOSMea_FunctionNet, 42
- SetSourceGate
 - CCMOSMea_FunctionNet, 42
- SetSpeedI
 - CRoboStatorDeviceNet, 391
- SetSpeedNativeI
 - CRoboStatorDeviceNet, 391
- SetSpeedNativeXY
 - CRoboStatorDeviceNet, 391
- SetSpeedNativeZ
 - CRoboStatorDeviceNet, 391
- SetSpeedXY
 - CRoboStatorDeviceNet, 391
- SetSpeedZ
 - CRoboStatorDeviceNet, 392
- SetStartTriggerSlope
 - CDigOutStimulatorFunctionNet, 60
- SetStgProgramInfo
 - CStg200xBasicNet, 443
- SetStimulusSites
 - CCMOSMea_FunctionNet, 42
- SetStopTriggerSlope
 - CDigOutStimulatorFunctionNet, 60
- SetStringFormat
 - CMcsUsbListEntryNet, 220
 - CMcsUsbListNet, 223
- SetSubChannel
 - CMcsBus_MotorControlNet, 135
- SetSwitches

- CSafeISDeviceNet, 395
- SetSyncoutMap
 - CStg200xBasicNet, 443
- SetTableName
 - CWarnerValveControllerDeviceNet, 549
- SetTablepointer
 - CRetinaLedDeviceNet, 345
- SetTableStep
 - CWarnerValveControllerDeviceNet, 549
- SetTableStepAll
 - CWarnerValveControllerDeviceNet, 549
- SetTestMode
 - CRFFFunctionNet, 350
- SetThermocoupleNanovoltPerKelvin
 - CFluidControlDeviceNet, 77
 - CTcxDeviceNet, 487
- SetThermoOffset
 - CMcsBus_TempSensorNet, 148
- SetTransformer
 - CMeFunctionNet, 280
- SetTrigger
 - CRetinaLedDeviceNet, 345
 - CWarnerValveControllerDeviceTesterFunctionNet, 557
- SetTriggerMaskValue
 - CMeaDeviceNet, 266
- SetTriggerPeriod
 - CMeaDeviceNet, 267
- SetTriggerSource
 - CStg200xBasicNet, 443, 444
- SetTriggerSyncDirection
 - CWarnerValveControllerDeviceTesterFunctionNet, 558
- SetUByteBuffer
 - CGenericDevelopDeviceNet, 94
- SetUClamp
 - CRoboDacqNet, 366
- SetUOffset
 - CRoboDacqNet, 366
- SetUIntBuffer
 - CGenericDevelopDeviceNet, 94
- SetupGroupDacqQueue
 - CMcsUsbDacqNet, 199
- SetupRetriggerMode
 - CStg200xDownloadBasicNet, 451
- SetupTrigger
 - CStg200xDownloadBasicNet, 451
 - CStimulusFunctionNet, 469
- SetupTriggerSingle
 - CStg200xDownloadBasicNet, 452
 - CStimulusFunctionNet, 469
- SetUseBubble
 - CPPS_FunctionNet, 335
- SetUserParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, 589
- SetUShortBuffer
 - CGenericDevelopDeviceNet, 95
- SetUOffset
 - CRoboDacqNet, 366
- SetValue
 - CGenericDevelopDeviceNet, 95
- SetValve
 - CFluidControlDeviceNet, 77
 - CRoboFluidDeviceNet, 384
- SetValve1
 - CRobo_FYIProgram_FunctionNet, 351
- SetValve2
 - CRobo_FYIProgram_FunctionNet, 352
- SetValveActive
 - CPPCFunctionNet, 330
 - CWarnerValveControllerDeviceNet, 549
- SetValveDigitalInInvert
 - CWarnerValveControllerDeviceNet, 550
- SetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, 550
- SetValveManualGroup
 - CWarnerValveControllerDeviceNet, 550
- SetValveManualState
 - CWarnerValveControllerDeviceNet, 550
- SetValveMode
 - CWarnerValveControllerDeviceNet, 551
- SetValves
 - CMcsBus_FYIExtensionNet, 118
- SetValvesActiveMap
 - CWarnerValveControllerDeviceNet, 551
- SetValvesManualStateMap
 - CWarnerValveControllerDeviceNet, 551
- SetValveTableEntry
 - CWarnerValveControllerDeviceNet, 551
- SetVelocityI
 - CRoboStatorDeviceNet, 392
- SetVelocityXY
 - CRoboStatorDeviceNet, 392
- SetVelocityZ
 - CRoboStatorDeviceNet, 392
- SetVMMaxNegativeCurrent
 - CMcsBus_VoltageModeNet, 151
- SetVMMaxNegativeCurrentEeprom
 - CMcsBus_VoltageModeNet, 151
- SetVMMaxNegativeVoltage
 - CMcsBus_VoltageModeNet, 151
- SetVMMaxNegativeVoltageEeprom
 - CMcsBus_VoltageModeNet, 151
- SetVMMaxPositiveCurrent
 - CMcsBus_VoltageModeNet, 152
- SetVMMaxPositiveCurrentEeprom
 - CMcsBus_VoltageModeNet, 152
- SetVMMaxPositiveVoltage
 - CMcsBus_VoltageModeNet, 152
- SetVMMaxPositiveVoltageEeprom
 - CMcsBus_VoltageModeNet, 152
- SetVMOutputOnOff
 - CMcsBus_VoltageModeNet, 152
- SetVMVoltage
 - CMcsBus_VoltageModeNet, 152

- SetVoltage12VLimit
 - CRoboDeviceNet, [376](#)
- SetVoltage5VLimit
 - CRoboDeviceNet, [376](#)
- SetVoltageAirvalveLimit
 - CRoboDeviceNet, [376](#)
- SetVoltageClampControllerParam_D
 - CWarnerUssingFunctionNet, [530](#)
- SetVoltageClampControllerParam_I
 - CWarnerUssingFunctionNet, [531](#)
- SetVoltageClampControllerParam_P
 - CWarnerUssingFunctionNet, [531](#)
- SetVoltageMode
 - CStg200xBasicNet, [444](#)
- SetVoltageRangeByIndex
 - CMcsUsbDacqNet, [199](#)
- SetVoltageRangeInMicroVolt
 - CMcsUsbDacqNet, [199](#)
- SetVoltageRs485ALimit
 - CRoboDeviceNet, [376](#)
- SetVoltageRs485BLimit
 - CRoboDeviceNet, [376](#)
- SetVoltageValvesLimit
 - CRoboDeviceNet, [376](#)
- SetWaveform
 - CChannelTestDeviceNet, [32](#)
 - CTEERFunctionNet, [496](#)
- SetWaveLengthInNanometer
 - CMultiwellOptoStimFunctionNet, [301](#)
- SetWorkingFrequency
 - CRFFFunctionNet, [350](#)
- SetWPADebugMode
 - CWClassicFunctionNet, [562](#)
- SetWPAType
 - CWClassicFunctionNet, [562](#)
- SetXGain
 - CRoboDacqNet, [366](#)
- Sideband
 - CStimulusFunctionNet::SidebandData, [591](#)
- SidebandData
 - CStimulusFunctionNet::SidebandData, [591](#)
- SineStart
 - COctoPotDeviceNet, [306](#)
- size
 - DigitalSource< digitalsourceenum >, [566](#)
- SN
 - HeadStageIDType, [577](#)
- Source
 - DigitalSource< digitalsourceenum >, [566](#)
- Start
 - CMeaCleanDeviceNet, [254](#)
 - CMeaCoatDeviceNet, [259](#)
 - CRobo_FYIProgram_FunctionNet, [352](#)
- StartDacq
 - CMcsUsbDacqNet, [199–201](#)
- StartInternalCalibration
 - CTEERFunctionNet, [496](#)
- StartLoop
 - CMcsUsbDacqNet, [201, 202](#)
- StartMCMovement
 - CMcsBus_MotorControlNet, [136](#)
- StartMeasurement
 - CMeaImpedanceDeviceNet, [275](#)
- StartPoll
 - CStimulusFunctionNet, [471](#)
 - CW2100_StimulatorFunctionNet, [512](#)
- StartSampling
 - CTEERFunctionNet, [496](#)
- StartSync
 - CMcsBus_SensorNet, [146](#)
- State
 - HeadStageIDTypeState, [579](#)
- Status
 - CUsbExceptionNet, [499](#)
- Status_AlreadyConfigured
 - CMcsUsbNet, [240](#)
- Status_BadStartFrame
 - CMcsUsbNet, [241](#)
- Status_Btstuff
 - CMcsUsbNet, [241](#)
- Status_BufferOverrun
 - CMcsUsbNet, [241](#)
- Status_BufferUnderrun
 - CMcsUsbNet, [241](#)
- Status_Canceled
 - CMcsUsbNet, [241](#)
- Status_Canceling
 - CMcsUsbNet, [241](#)
- Status_ConnectedPipes
 - CMcsUsbNet, [241](#)
- Status_ControlNotOwned
 - CMcsUsbNet, [241](#)
- Status_Crc
 - CMcsUsbNet, [241](#)
- Status_DataOverrun
 - CMcsUsbNet, [241](#)
- Status_DataToggleMismatch
 - CMcsUsbNet, [241](#)
- Status_DataUnderrun
 - CMcsUsbNet, [242](#)
- Status_DeviceLocked
 - CMcsUsbNet, [242](#)
- Status_DeviceNotFound
 - CMcsUsbNet, [242](#)
- Status_DeviceRemoved
 - CMcsUsbNet, [242](#)
- Status_DevNotResponding
 - CMcsUsbNet, [242](#)
- Status_EndpointHalted
 - CMcsUsbNet, [242](#)
- Status_ErrorBusy
 - CMcsUsbNet, [242](#)
- Status_ErrorShortTransfer
 - CMcsUsbNet, [242](#)
- Status_Fifo
 - CMcsUsbNet, [242](#)

- Status_FrameControlOwned
 - CMcsUsbNet, [242](#)
- Status_InternalHcError
 - CMcsUsbNet, [242](#)
- Status_InvalidParameter
 - CMcsUsbNet, [243](#)
- Status_InvalidPipeHandle
 - CMcsUsbNet, [243](#)
- Status_InvalidUrbFunction
 - CMcsUsbNet, [243](#)
- Status_IoPending
 - CMcsUsbNet, [243](#)
- Status_IoTimeout
 - CMcsUsbNet, [243](#)
- Status_IsochRequestFailed
 - CMcsUsbNet, [243](#)
- Status_LastUsbErrorMismatch
 - CMcsUsbNet, [243](#)
- Status_NoBandwidth
 - CMcsUsbNet, [243](#)
- Status_NoMemory
 - CMcsUsbNet, [243](#)
- Status_NoSuchDevice
 - CMcsUsbNet, [243](#)
- Status_NotAccessed
 - CMcsUsbNet, [243](#)
- Status_NotSupported
 - CMcsUsbNet, [244](#)
- Status_PidCheckFailure
 - CMcsUsbNet, [244](#)
- Status_PipeNotLinked
 - CMcsUsbNet, [244](#)
- Status_RequestFailed
 - CMcsUsbNet, [244](#)
- Status_RequestMutexFailed
 - CMcsUsbNet, [244](#)
- Status_RequestMutexTimeout
 - CMcsUsbNet, [244](#)
- Status_Stall
 - CMcsUsbNet, [244](#)
- Status_Unconfigured
 - CMcsUsbNet, [244](#)
- Status_UnexpectedPid
 - CMcsUsbNet, [244](#)
- Stg200xPollStatusEvent
 - CStg200xDownloadNet, [459](#)
- StgStatusNet, [591](#)
 - FromIntPtr, [592](#)
 - FromPtr, [592](#)
 - ListOfChangedTriggers, [592](#)
 - TiggerStatus, [592](#)
- StillConnected
 - CRadioControlledDevicesNet, [342](#)
- Stimulator
 - CW2100_FunctionNet, [507](#)
- Stimulus
 - CCMOSMeaDeviceNet, [45](#)
 - CStg200xDownloadBasicNet, [453](#)
- StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [592](#)
- StimulusFunction
 - CLIH3DeviceNet, [114](#)
- StimulusParameters
 - HeadStageIDType, [577](#)
- Stop
 - CMeaCleanDeviceNet, [254](#)
 - CMeaCoatDeviceNet, [260](#)
- StopDacq
 - CMcsUsbDacqNet, [203](#)
- StopLoop
 - CMcsUsbDacqNet, [204](#)
- StopMCMovement
 - CMcsBus_MotorControlNet, [136](#)
- StopMovement
 - CRoboDeviceNet, [376](#)
- StopMovementI
 - CRoboStatorDeviceNet, [392](#)
- StopMovementXY
 - CRoboStatorDeviceNet, [392](#)
- StopMovementZ
 - CRoboStatorDeviceNet, [392](#)
- StopPlateClamp
 - CMultiwellDeviceNet, [296](#)
- StopPoll
 - CStimulusFunctionNet, [471](#)
 - CW2100_StimulatorFunctionNet, [513](#)
- StopSampling
 - CTEERFunctionNet, [496](#)
- StopTable
 - CRoboDacqNet, [366](#)
- StoreValveTable
 - CWarnerValveControllerDeviceNet, [552](#)
- SwitchOnOff
 - CPositionIIDeviceNet, [319](#)
- SYNC_BIT0
 - CW2100_StimulatorFunctionNet, [513](#)
- SYNC_BIT1
 - CW2100_StimulatorFunctionNet, [513](#)
- Table_Wait
 - CRoboDacqNet, [367](#)
- TableDefBegin
 - CRoboDacqNet, [367](#)
- TableDefEnd
 - CRoboDacqNet, [367](#)
- TableEntryChangedEvent
 - CWarnerValveControllerDeviceNet, [555](#)
- TactSwitchGetState
 - CMcsBus_SensorNet, [146](#)
- TactSwitchSetDisplay
 - CMcsBus_SensorNet, [146](#)
- TEERFunctionNet
 - CTEERMachineDeviceNet, [497](#)
- ThrowCUsbExceptionNetOnError
 - CMcsUsbFunctionNet, [216](#)
 - CMcsUsbNet, [238](#)

- TiggerStatus
 - StgStatusNet, [592](#)
- TimeResolutionInNanoSeconds
 - W2100_StimulusParametersNet, [595](#)
- ToString
 - CFilterPropertyNet, [71](#)
 - CMcsUsbListEntryNet, [220](#)
 - HeadStageIDType, [576](#)
 - HeadstageIDTypeObject, [578](#)
- TriggerStatus
 - CMcsUsbDeviceStatePushFunctionNet, [206](#)
 - CMcsUsbDeviceStatePushNet, [207](#)
- TxnGetSerialNumber
 - CMcsUsbNet, [238](#)
- TxnSetSerialNumber
 - CMcsUsbNet, [238](#)
- TxnTestMemoryReadAndCheck
 - CMcsUsbNet, [238](#)
- TxnTestMemoryWrite
 - CMcsUsbNet, [239](#)
- Type
 - HeadStageIDType, [577](#)
- TypeValue
 - HeadStageIDType, [577](#)
- UIntA1
 - CFilterCoefficientsNet, [66](#)
- UIntA2
 - CFilterCoefficientsNet, [66](#)
- UIntB0
 - CFilterCoefficientsNet, [66](#)
- UIntB1
 - CFilterCoefficientsNet, [66](#)
- UIntB2
 - CFilterCoefficientsNet, [66](#)
- Unknown
 - HeadStageIDType, [576](#)
- unknown
 - Mcs::Usb, [25](#)
- UnlockPlateClamp
 - CMultiwellDeviceNet, [296](#)
- UnrolledAmplitude
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- UnrolledDuration
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- UnrolledSync
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [593](#)
- UpdateChannelBlock
 - CCMOSMeaDeviceNet, [45](#)
- UpdateDisplay
 - CRoboDacqNet, [367](#)
- UpdateFirmware
 - CMcsUsbFactoryNet, [213](#), [214](#)
- UpdateTransistorVoltages
 - CCMOSMea_FunctionNet, [42](#)
- USB
 - FirmwareDestinationNames, [575](#)
 - usbSetupPacket_t, [593](#)
 - bmRequestType, [594](#)
 - bRequest, [594](#)
 - wIndex, [594](#)
 - wLength, [594](#)
 - wValue, [594](#)
 - UserDefinedName
 - HeadStageIDType, [577](#)
 - Valid
 - HeadStageIDType, [577](#)
 - ValidKey
 - CMcsUsbNet, [239](#)
 - Voltage
 - BatteryState, [28](#)
 - VoltageRangeDisplayStringMilliVolt
 - CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, [499](#)
 - VoltageRangeInMicroVolt
 - CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, [499](#)
 - W2100_StimulusParametersNet, [595](#)
 - VoltageResolutionInMicroVolt
 - W2100_StimulusParametersNet, [595](#)
 - VoltageString
 - BatteryState, [28](#)
 - VOPSTimerSetResetTimes
 - CCMOSMea_FunctionNet, [42](#), [43](#)
 - W16lsW14
 - HeadStageIDType, [577](#)
 - W2100_FunctionNet
 - CMeaDeviceNet, [268](#)
 - W2100_StimulusParametersNet, [594](#)
 - CurrentRangeInNanoAmp, [594](#)
 - CurrentResolutionInNanoAmp, [594](#)
 - DACResolution, [594](#)
 - TimeResolutionInNanoSeconds, [595](#)
 - VoltageRangeInMicroVolt, [595](#)
 - VoltageResolutionInMicroVolt, [595](#)
 - WaitForAllChambers
 - CWarnerUssingFunctionNet, [531](#)
 - WaitForChamber
 - CWarnerUssingFunctionNet, [531](#)
 - WaitForUssingFunction
 - CWarnerUssingDeviceNet, [515](#)
 - WClassicFunctionNet
 - CMeaDeviceNet, [268](#)
 - wIndex
 - usbSetupPacket_t, [594](#)
 - wLength
 - usbSetupPacket_t, [594](#)
 - WPAError_ScanningIsPending
 - CMcsUsbNet, [244](#)
 - Write
 - CExternDTesterDeviceNet, [63](#)
 - Write2
 - CExternDTesterDeviceNet, [63](#)

WriteEepromRegisterPreconfig
 CMcsUsbNet, [239](#)
WritePipe
 CGenericDevelopDeviceNet, [96](#)
WriteRegister
 CMcsUsbNet, [239](#), [240](#)
WriteRegister32
 CMcsUsbNet, [240](#)
WriteRegisterArray
 CMcsUsbNet, [240](#)
WriteRegisterTimeSlot
 CMcsUsbNet, [240](#)
WriteRegisterValue
 CMcsUsbNet, [240](#)
wValue
 usbSetupPacket_t, [594](#)