

multichannel* systems

McsUsbNet.dll

Version 5.1.25

Multi Channel Systems MCS GmbH
Aspenhastrasse 21
72770 Reutlingen
Germany
Fon +49-71 21-90 92 5 - 0
Fax +49-71 21-90 92 5 -11
info@multichannelsystems.com
www.multichannelsystems.com

Generated by Doxygen 1.9.1

1 McsUsbNet.dll for MCS USB devices	1
1.1 Introduction	1
1.2 System requirements	1
1.3 Connecting to an MCS device	2
2 Device Classes	2
2.1 The MCS FluidControl Device	2
2.1.1 Introduction	2
2.1.2 Access to the FluidControl device	2
2.2 MCS-USB-Sw2to64 device	3
3 Function Classes	3
4 Data ACQuisition (DACQ) Devices	4
5 The MCS Robo Device	5
5.1 Introduction	5
6 STG200x & STG400x STimulus Generator	5
6.1 Introduction	5
6.2 Download mode	6
6.2.1 Memory Layout and Trigger Setup	6
6.3 Streaming mode	8
6.3.1 Memory Layout and Trigger Setup	9
7 Namespace Index	11
7.1 Namespace List	11
8 Hierarchical Index	11
8.1 Class Hierarchy	11
9 Class Index	16
9.1 Class List	16
10 Namespace Documentation	22
10.1 Mcs Namespace Reference	22
10.2 Mcs::Usb Namespace Reference	22
10.2.1 Enumeration Type Documentation	50
10.2.2 Function Documentation	90
11 Class Documentation	92
11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference	92
11.1.1 Member Data Documentation	92
11.2 BatteryState Class Reference	92
11.2.1 Property Documentation	92
11.3 BesselFilterHighPassNet Class Reference	93

11.3.1 Constructor & Destructor Documentation	93
11.4 BesselFilterLowPassNet Class Reference	93
11.4.1 Constructor & Destructor Documentation	94
11.5 ButterworthFilterHighPassNet Class Reference	94
11.5.1 Constructor & Destructor Documentation	94
11.6 ButterworthFilterLowPassNet Class Reference	95
11.6.1 Constructor & Destructor Documentation	95
11.7 CChannelTestDeviceNet Class Reference	95
11.7.1 Constructor & Destructor Documentation	96
11.7.2 Member Function Documentation	96
11.8 CCMOSMea_FunctionNet Class Reference	96
11.8.1 Constructor & Destructor Documentation	98
11.8.2 Member Function Documentation	99
11.9 CCMOSMeaDeviceNet Class Reference	107
11.9.1 Constructor & Destructor Documentation	108
11.9.2 Member Function Documentation	108
11.9.3 Property Documentation	110
11.10 CCreateFilterNet Class Reference	110
11.10.1 Constructor & Destructor Documentation	111
11.10.2 Member Function Documentation	111
11.10.3 Property Documentation	112
11.11 CDacCalibrationFunctionNet Class Reference	112
11.11.1 Detailed Description	113
11.11.2 Constructor & Destructor Documentation	113
11.11.3 Member Function Documentation	113
11.12 CDacqGroupChannelGenericSelectionNet Class Reference	114
11.12.1 Constructor & Destructor Documentation	115
11.13 CDacqGroupChannelSelectionNet Class Reference	115
11.13.1 Constructor & Destructor Documentation	115
11.14 CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference	115
11.14.1 Constructor & Destructor Documentation	116
11.14.2 Member Function Documentation	116
11.15 CDeviceGroupChannelInfoGenericNet Class Reference	118
11.15.1 Constructor & Destructor Documentation	118
11.16 CDeviceGroupChannelInfoMEA2100_256Net Class Reference	119
11.16.1 Constructor & Destructor Documentation	119
11.17 CDeviceGroupChannelInfoNet Class Reference	119
11.17.1 Constructor & Destructor Documentation	120
11.18 CDeviceGroupChannelInfoSCUNet Class Reference	120
11.18.1 Constructor & Destructor Documentation	120
11.19 CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet > Class Template Reference	120

11.19.1 Constructor & Destructor Documentation	121
11.19.2 Member Data Documentation	121
11.20 CDeviceGroupChannelInfoW2100Net Class Reference	121
11.20.1 Constructor & Destructor Documentation	122
11.21 CDigOutStimulatorFunctionNet Class Reference	122
11.21.1 Detailed Description	123
11.21.2 Constructor & Destructor Documentation	123
11.21.3 Member Function Documentation	123
11.22 CEncapsulatorDeviceNet Class Reference	126
11.22.1 Detailed Description	126
11.22.2 Constructor & Destructor Documentation	127
11.22.3 Member Function Documentation	127
11.23 CExternDTesterDeviceNet Class Reference	127
11.23.1 Detailed Description	127
11.23.2 Constructor & Destructor Documentation	128
11.23.3 Member Function Documentation	128
11.24 CFilterCoefficientsNet Class Reference	129
11.24.1 Constructor & Destructor Documentation	129
11.24.2 Member Function Documentation	130
11.24.3 Property Documentation	131
11.25 CFilterConfigurationNet Class Reference	131
11.25.1 Constructor & Destructor Documentation	131
11.25.2 Member Function Documentation	132
11.26 CFilterConfigurationRegisterNet Class Reference	133
11.26.1 Constructor & Destructor Documentation	133
11.26.2 Member Function Documentation	133
11.27 CFilterPropertyNet Class Reference	135
11.27.1 Constructor & Destructor Documentation	135
11.27.2 Member Function Documentation	135
11.27.3 Property Documentation	135
11.28 CFluidControlDeviceNet Class Reference	136
11.28.1 Detailed Description	137
11.28.2 Constructor & Destructor Documentation	138
11.28.3 Member Function Documentation	138
11.28.4 Property Documentation	143
11.29 CFYIDeviceNet Class Reference	143
11.29.1 Detailed Description	144
11.29.2 Constructor & Destructor Documentation	144
11.29.3 Property Documentation	144
11.30 CGenericDevelopDeviceNet Class Reference	145
11.30.1 Detailed Description	152
11.30.2 Constructor & Destructor Documentation	152

11.30.3 Member Function Documentation	152
11.31 CGilsonDeviceNet Class Reference	163
11.31.1 Detailed Description	163
11.31.2 Constructor & Destructor Documentation	163
11.31.3 Member Function Documentation	164
11.31.4 Member Data Documentation	164
11.32 CGrapheneASICDeviceNet Class Reference	165
11.32.1 Constructor & Destructor Documentation	165
11.32.2 Member Function Documentation	165
11.33 CGrapheneFunctionNet Class Reference	166
11.33.1 Detailed Description	168
11.33.2 Constructor & Destructor Documentation	168
11.33.3 Member Function Documentation	168
11.34 CHiClampDeviceNet Class Reference	177
11.34.1 Detailed Description	177
11.34.2 Constructor & Destructor Documentation	177
11.34.3 Property Documentation	178
11.35 CHLADacqNet Class Reference	178
11.35.1 Constructor & Destructor Documentation	178
11.36 CHLADeviceNet Class Reference	178
11.36.1 Detailed Description	179
11.36.2 Constructor & Destructor Documentation	179
11.36.3 Property Documentation	179
11.37 CMcsUsbDacqNet::CHWInfo Class Reference	179
11.37.1 Detailed Description	180
11.37.2 Constructor & Destructor Documentation	180
11.37.3 Member Function Documentation	180
11.38 CIntanMea_FunctionNet Class Reference	182
11.38.1 Constructor & Destructor Documentation	183
11.38.2 Member Function Documentation	183
11.39 CInterfaceboard2FunctionNet Class Reference	184
11.39.1 Detailed Description	185
11.39.2 Constructor & Destructor Documentation	185
11.39.3 Member Function Documentation	185
11.40 CInterfaceboardFunctionNet Class Reference	186
11.40.1 Detailed Description	187
11.40.2 Constructor & Destructor Documentation	187
11.40.3 Member Function Documentation	187
11.41 CLIH3DeviceNet Class Reference	188
11.41.1 Detailed Description	190
11.41.2 Constructor & Destructor Documentation	190
11.41.3 Member Function Documentation	190

11.41.4 Property Documentation	197
11.42 CMcsBus_AxisParametersNet Class Reference	197
11.42.1 Constructor & Destructor Documentation	198
11.42.2 Member Function Documentation	198
11.43 CMcsBus_ExtensionNet Class Reference	199
11.43.1 Constructor & Destructor Documentation	199
11.43.2 Member Function Documentation	200
11.44 CMcsBus_FYIExtensionNet Class Reference	200
11.44.1 Constructor & Destructor Documentation	200
11.44.2 Member Function Documentation	201
11.45 CMcsBus_MotorControlNet Class Reference	202
11.45.1 Constructor & Destructor Documentation	205
11.45.2 Member Function Documentation	205
11.46 CMcsBus_SensorNet Class Reference	219
11.46.1 Constructor & Destructor Documentation	221
11.46.2 Member Function Documentation	221
11.47 CMcsBus_TempSensorNet Class Reference	229
11.47.1 Constructor & Destructor Documentation	230
11.47.2 Member Function Documentation	230
11.48 CMcsBus_VoltageModeNet Class Reference	231
11.48.1 Constructor & Destructor Documentation	232
11.48.2 Member Function Documentation	233
11.49 CMcsBusNet Class Reference	236
11.49.1 Constructor & Destructor Documentation	236
11.49.2 Member Function Documentation	237
11.50 CMcsUsbDacqNet Class Reference	240
11.50.1 Detailed Description	246
11.50.2 Constructor & Destructor Documentation	246
11.50.3 Member Function Documentation	246
11.50.4 Member Data Documentation	286
11.50.5 Property Documentation	287
11.50.6 Event Documentation	287
11.51 CMcsUsbDeviceStatePushFunctionNet Class Reference	287
11.51.1 Constructor & Destructor Documentation	288
11.51.2 Member Function Documentation	288
11.51.3 Event Documentation	288
11.52 CMcsUsbDeviceStatePushNet Class Reference	288
11.52.1 Constructor & Destructor Documentation	289
11.52.2 Member Function Documentation	289
11.52.3 Event Documentation	289
11.53 CMcsUsbFactoryNet Class Reference	289
11.53.1 Constructor & Destructor Documentation	291

11.53.2 Member Function Documentation	291
11.53.3 Member Data Documentation	297
11.54 CMcsUsbFunctionNet Class Reference	298
11.54.1 Constructor & Destructor Documentation	298
11.54.2 Member Function Documentation	299
11.54.3 Member Data Documentation	299
11.55 CMcsUsbFunctionPointerContainer Class Reference	299
11.56 CMcsUsbListEntryNet Class Reference	299
11.56.1 Detailed Description	300
11.56.2 Constructor & Destructor Documentation	300
11.56.3 Member Function Documentation	301
11.56.4 Property Documentation	304
11.57 CMcsUsbListNet Class Reference	305
11.57.1 Detailed Description	305
11.57.2 Constructor & Destructor Documentation	305
11.57.3 Member Function Documentation	306
11.57.4 Property Documentation	307
11.57.5 Event Documentation	307
11.58 CMcsUsbNet Class Reference	308
11.58.1 Detailed Description	312
11.58.2 Constructor & Destructor Documentation	312
11.58.3 Member Function Documentation	312
11.58.4 Member Data Documentation	324
11.58.5 Property Documentation	328
11.59 CMcsUsbPointerContainer Class Reference	329
11.60 CMEA2100_256DacqGroupChannelSelectionNet Class Reference	329
11.60.1 Constructor & Destructor Documentation	329
11.61 CMEA2100x256FunctionNet Class Reference	329
11.61.1 Detailed Description	330
11.61.2 Constructor & Destructor Documentation	330
11.61.3 Member Function Documentation	330
11.62 CMeaAudioFunctionNet Class Reference	331
11.62.1 Constructor & Destructor Documentation	331
11.62.2 Member Function Documentation	332
11.63 CMeaCleanDeviceNet Class Reference	334
11.63.1 Detailed Description	335
11.63.2 Constructor & Destructor Documentation	335
11.63.3 Member Function Documentation	335
11.64 CMeaCoatDeviceNet Class Reference	338
11.64.1 Detailed Description	339
11.64.2 Constructor & Destructor Documentation	339
11.64.3 Member Function Documentation	339

11.65 CMeaDeviceNet Class Reference	343
11.65.1 Detailed Description	344
11.65.2 Constructor & Destructor Documentation	345
11.65.3 Member Function Documentation	345
11.65.4 Property Documentation	350
11.66 CMeaDigitalDataFunctionNet Class Reference	351
11.66.1 Constructor & Destructor Documentation	352
11.66.2 Member Function Documentation	352
11.67 CMeaFeedbackFunctionNet Class Reference	353
11.67.1 Constructor & Destructor Documentation	354
11.67.2 Member Function Documentation	354
11.68 CMeaImpedanceDeviceNet Class Reference	357
11.68.1 Constructor & Destructor Documentation	358
11.68.2 Member Function Documentation	358
11.69 CMeasureTableDeviceNet Class Reference	359
11.69.1 Detailed Description	359
11.69.2 Constructor & Destructor Documentation	359
11.69.3 Property Documentation	359
11.70 CMeaSwitchDeviceNet Class Reference	360
11.70.1 Detailed Description	360
11.70.2 Constructor & Destructor Documentation	360
11.70.3 Member Function Documentation	361
11.71 CMeaUSBDeviceNet Class Reference	361
11.71.1 Detailed Description	362
11.71.2 Constructor & Destructor Documentation	362
11.72 CMeFunctionNet Class Reference	363
11.72.1 Detailed Description	363
11.72.2 Constructor & Destructor Documentation	363
11.72.3 Member Function Documentation	364
11.73 CMultiBatteryChargerDeviceNet Class Reference	364
11.73.1 Detailed Description	365
11.73.2 Constructor & Destructor Documentation	365
11.73.3 Member Function Documentation	365
11.74 CMultiwellCallbackFunctionNet Class Reference	372
11.74.1 Detailed Description	373
11.74.2 Constructor & Destructor Documentation	373
11.74.3 Member Function Documentation	373
11.74.4 Event Documentation	374
11.75 CMultiwellDeviceNet Class Reference	374
11.75.1 Detailed Description	375
11.75.2 Constructor & Destructor Documentation	375
11.75.3 Member Function Documentation	376

11.76 CMultiwellOptoStimFunctionNet Class Reference	382
11.76.1 Detailed Description	382
11.76.2 Constructor & Destructor Documentation	382
11.76.3 Member Function Documentation	383
11.77 CNF_GenDeviceNet Class Reference	386
11.77.1 Constructor & Destructor Documentation	387
11.77.2 Member Function Documentation	387
11.78 COctoPotDeviceNet Class Reference	387
11.78.1 Constructor & Destructor Documentation	388
11.78.2 Member Function Documentation	388
11.79 COkuvisionStimulatorDeviceNet Class Reference	391
11.79.1 Constructor & Destructor Documentation	392
11.79.2 Member Function Documentation	392
11.80 CPatchServerDeviceNet Class Reference	395
11.80.1 Detailed Description	395
11.80.2 Constructor & Destructor Documentation	395
11.80.3 Property Documentation	395
11.81 CPathIdentDeviceNet Class Reference	396
11.81.1 Constructor & Destructor Documentation	396
11.81.2 Member Function Documentation	396
11.82 CPedoterDeviceNet Class Reference	397
11.82.1 Detailed Description	397
11.82.2 Constructor & Destructor Documentation	397
11.82.3 Member Function Documentation	397
11.83 CPeristalticPumpDeviceNet Class Reference	398
11.83.1 Detailed Description	398
11.83.2 Constructor & Destructor Documentation	399
11.83.3 Property Documentation	399
11.84 CPgaDeviceNet Class Reference	399
11.84.1 Constructor & Destructor Documentation	400
11.84.2 Member Function Documentation	400
11.85 CPositionIIDeviceNet Class Reference	401
11.85.1 Detailed Description	403
11.85.2 Constructor & Destructor Documentation	403
11.85.3 Member Function Documentation	403
11.85.4 Property Documentation	410
11.86 CPositionImpDeviceNet Class Reference	410
11.86.1 Detailed Description	410
11.86.2 Constructor & Destructor Documentation	411
11.86.3 Member Function Documentation	411
11.87 CPPCDeviceNet Class Reference	413
11.87.1 Constructor & Destructor Documentation	413

11.87.2 Property Documentation	414
11.88 CPPCFunctionNet Class Reference	414
11.88.1 Detailed Description	415
11.88.2 Constructor & Destructor Documentation	415
11.88.3 Member Function Documentation	416
11.89 CPPS_DeviceNet Class Reference	423
11.89.1 Constructor & Destructor Documentation	423
11.89.2 Property Documentation	423
11.90 CPPS_FunctionNet Class Reference	424
11.90.1 Constructor & Destructor Documentation	424
11.90.2 Member Function Documentation	425
11.91 CPPSDeviceNet Class Reference	428
11.91.1 Detailed Description	428
11.91.2 Constructor & Destructor Documentation	429
11.92 CProgramPressureCurveNet Class Reference	429
11.92.1 Detailed Description	429
11.92.2 Constructor & Destructor Documentation	429
11.92.3 Member Function Documentation	430
11.93 CPulseGeneratorFunctionNet Class Reference	430
11.93.1 Detailed Description	431
11.93.2 Constructor & Destructor Documentation	431
11.93.3 Member Function Documentation	432
11.94 CRadioControlledDevicesNet Class Reference	433
11.94.1 Constructor & Destructor Documentation	434
11.94.2 Member Function Documentation	434
11.95 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference	435
11.95.1 Constructor & Destructor Documentation	435
11.95.2 Member Function Documentation	435
11.95.3 Member Data Documentation	436
11.96 CRetinaLedDeviceNet Class Reference	436
11.96.1 Constructor & Destructor Documentation	437
11.96.2 Member Function Documentation	437
11.97 CRFFFunctionNet Class Reference	438
11.97.1 Detailed Description	439
11.97.2 Constructor & Destructor Documentation	439
11.97.3 Member Function Documentation	439
11.98 CRobo_FYIPProgram_FunctionNet Class Reference	443
11.98.1 Constructor & Destructor Documentation	443
11.98.2 Member Function Documentation	443
11.99 CRobo_FYITemp_FunctionNet Class Reference	444
11.99.1 Constructor & Destructor Documentation	445
11.99.2 Member Function Documentation	445

11.100 CRoboDacqNet Class Reference	446
11.100.1 Constructor & Destructor Documentation	449
11.100.2 Member Function Documentation	449
11.100.3 Member Data Documentation	460
11.101 CRoboDeviceNet Class Reference	461
11.101.1 Detailed Description	464
11.101.2 Constructor & Destructor Documentation	464
11.101.3 Member Function Documentation	464
11.101.4 Member Data Documentation	472
11.101.5 Property Documentation	476
11.101.6 Event Documentation	477
11.102 CRoboFluidDeviceNet Class Reference	477
11.102.1 Constructor & Destructor Documentation	478
11.102.2 Member Function Documentation	478
11.102.3 Member Data Documentation	480
11.102.4 Property Documentation	480
11.103 CRoboInjectDeviceNet Class Reference	480
11.103.1 Detailed Description	481
11.103.2 Constructor & Destructor Documentation	481
11.104 CRoboocyte2DeviceNet Class Reference	481
11.104.1 Detailed Description	481
11.104.2 Constructor & Destructor Documentation	482
11.104.3 Member Function Documentation	482
11.105 CRoboStatorDeviceNet Class Reference	482
11.105.1 Constructor & Destructor Documentation	484
11.105.2 Member Function Documentation	484
11.105.3 Property Documentation	488
11.106 CSafeISDeviceNet Class Reference	488
11.106.1 Detailed Description	489
11.106.2 Constructor & Destructor Documentation	489
11.106.3 Member Function Documentation	490
11.106.4 Property Documentation	491
11.107 CSCUDacqGroupChannelSelectionNet Class Reference	491
11.107.1 Constructor & Destructor Documentation	492
11.108 CSCUFunctionNet Class Reference	492
11.108.1 Detailed Description	494
11.108.2 Constructor & Destructor Documentation	494
11.108.3 Member Function Documentation	495
11.108.4 Event Documentation	506
11.109 CSerialPortNet Class Reference	507
11.109.1 Constructor & Destructor Documentation	507
11.109.2 Member Function Documentation	507

11.110 CStg200xBasicNet Class Reference	508
11.110.1 Detailed Description	512
11.110.2 Constructor & Destructor Documentation	513
11.110.3 Member Function Documentation	513
11.111 CStg200xDownloadBasicNet Class Reference	542
11.111.1 Detailed Description	543
11.111.2 Member Function Documentation	543
11.111.3 Property Documentation	549
11.112 CStg200xDownloadNet Class Reference	549
11.112.1 Detailed Description	550
11.112.2 Constructor & Destructor Documentation	550
11.112.3 Member Function Documentation	551
11.112.4 Event Documentation	555
11.113 CStimulusFunctionNet Class Reference	556
11.113.1 Constructor & Destructor Documentation	557
11.113.2 Member Function Documentation	558
11.113.3 Event Documentation	567
11.114 CSw2to64DeviceNet Class Reference	567
11.114.1 Detailed Description	567
11.114.2 Constructor & Destructor Documentation	568
11.114.3 Member Function Documentation	568
11.115 CTcxDeviceNet Class Reference	569
11.115.1 Detailed Description	571
11.115.2 Constructor & Destructor Documentation	571
11.115.3 Member Function Documentation	572
11.116 CTEERFunctionNet Class Reference	582
11.116.1 Detailed Description	584
11.116.2 Constructor & Destructor Documentation	584
11.116.3 Member Function Documentation	585
11.117 CTEERMachineDeviceNet Class Reference	593
11.117.1 Constructor & Destructor Documentation	593
11.117.2 Property Documentation	594
11.118 CUsbDeviceConfigurationFunctionNet Class Reference	594
11.118.1 Detailed Description	594
11.118.2 Constructor & Destructor Documentation	594
11.118.3 Member Function Documentation	595
11.119 CUsbExceptionNet Class Reference	595
11.119.1 Detailed Description	596
11.119.2 Constructor & Destructor Documentation	596
11.119.3 Property Documentation	596
11.120 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference	596
11.120.1 Constructor & Destructor Documentation	597

11.120.2 Member Data Documentation	597
11.121 CW2100_FunctionNet Class Reference	597
11.121.1 Constructor & Destructor Documentation	599
11.121.2 Member Function Documentation	599
11.121.3 Property Documentation	604
11.122 CW2100_StimulatorFunctionNet Class Reference	605
11.122.1 Constructor & Destructor Documentation	606
11.122.2 Member Function Documentation	606
11.122.3 Member Data Documentation	611
11.122.4 Event Documentation	611
11.123 CW2100DacqGroupChannelSelectionNet Class Reference	611
11.123.1 Constructor & Destructor Documentation	611
11.124 CWarnerUssingDeviceNet Class Reference	612
11.124.1 Detailed Description	612
11.124.2 Constructor & Destructor Documentation	612
11.124.3 Property Documentation	613
11.125 CWarnerUssingFunctionNet Class Reference	613
11.125.1 Detailed Description	615
11.125.2 Constructor & Destructor Documentation	615
11.125.3 Member Function Documentation	615
11.126 CWarnerValveControllerDeviceNet Class Reference	629
11.126.1 Detailed Description	633
11.126.2 Constructor & Destructor Documentation	634
11.126.3 Member Function Documentation	634
11.126.4 Event Documentation	651
11.127 CWarnerValveControllerDeviceTesterFunctionNet Class Reference	654
11.127.1 Detailed Description	655
11.127.2 Constructor & Destructor Documentation	655
11.127.3 Member Function Documentation	656
11.128 CWClassicFunctionNet Class Reference	657
11.128.1 Constructor & Destructor Documentation	658
11.128.2 Member Function Documentation	659
11.129 CWirelessBaseFunctionNet Class Reference	662
11.129.1 Constructor & Destructor Documentation	663
11.129.2 Member Function Documentation	663
11.130 DeviceIdNet Struct Reference	663
11.130.1 Detailed Description	663
11.130.2 Constructor & Destructor Documentation	664
11.130.3 Member Function Documentation	664
11.130.4 Member Data Documentation	664
11.131 DigitalSource< digitalsourceenum > Class Template Reference	665
11.131.1 Constructor & Destructor Documentation	665

11.131.2 Member Function Documentation	665
11.131.3 Property Documentation	666
11.132 DigitalSourceGeneral Class Reference	666
11.132.1 Constructor & Destructor Documentation	666
11.132.2 Member Function Documentation	666
11.132.3 Property Documentation	667
11.133 DriverVersionNet Class Reference	667
11.133.1 Detailed Description	668
11.133.2 Constructor & Destructor Documentation	668
11.133.3 Member Function Documentation	668
11.134 FirmwareDestinationNames Class Reference	673
11.134.1 Member Data Documentation	674
11.135 HeadStageIDType Class Reference	676
11.135.1 Member Enumeration Documentation	677
11.135.2 Constructor & Destructor Documentation	677
11.135.3 Member Function Documentation	678
11.135.4 Property Documentation	678
11.136 HeadstageIDTypeObject Class Reference	679
11.136.1 Constructor & Destructor Documentation	680
11.136.2 Member Function Documentation	680
11.136.3 Member Data Documentation	680
11.136.4 Property Documentation	681
11.137 HeadStageIDTypeState Class Reference	681
11.137.1 Property Documentation	681
11.138 mkfilterNet Class Reference	682
11.138.1 Member Function Documentation	682
11.139 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference	685
11.139.1 Member Function Documentation	686
11.140 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference	692
11.140.1 Member Function Documentation	692
11.141 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference	692
11.141.1 Constructor & Destructor Documentation	692
11.141.2 Member Function Documentation	692
11.141.3 Member Data Documentation	693
11.142 CMeaAudioFunctionNet::s_setaudionet Struct Reference	693
11.142.1 Member Data Documentation	693
11.143 CStimulusFunctionNet::SidebandData Class Reference	694
11.143.1 Constructor & Destructor Documentation	694
11.143.2 Property Documentation	694
11.144 StgStatusNet Class Reference	695
11.144.1 Member Function Documentation	695
11.144.2 Member Data Documentation	695

11.145 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference	695
11.145.1 Constructor & Destructor Documentation	696
11.145.2 Property Documentation	696
11.146 usbSetupPacket_t Class Reference	697
11.146.1 Member Data Documentation	697
11.147 W2100_StimulusParametersNet Struct Reference	697
11.147.1 Member Data Documentation	698
Index	699

1 McsUsbNet.dll for MCS USB devices

1.1 Introduction

This DLL provides the .NET interface to MCS devices

The most important options are accessing our stimulator and data acquisition devices:

- [STG200x & STG400x STimulus Generator](#)
- [Data ACQuisition \(DACQ\) Devices](#)

See here for a list of our other devices: [Device Classes](#).

And here for a list of function classes addressing groups of features that might be shared between different devices: [Function Classes](#).

1.2 System requirements

The DLL can be used with any .NET compatible language.

The DLL needs the **.NET Framework 4.7.2**.

It requires the **Microsoft Visual C++ Redistributable for Visual Studio 2019** to be installed.

It also requires the **USB driver** to be installed.

The simplest way to achieve this is to install the latest **Multi Channel Experimenter** setup (will install 64bit redistributable).

All examples assume that the [Mcs.Usb](#) namespace is loaded:

```
using namespace Mcs.Usb;
```

Include the file McsUsbNet.dll into the references of your project.

1.3 Connecting to an MCS device

A connection to a DAQ device is established by [Mcs.Usb.CMcsUsbNet.Connect](#). When this function is called without argument, the first DAQ device found on the USB bus is used:

```
CMcsUsbNet device = new CMcsUsbNet();  
device.Connect();
```

When more than one DAQ device of the specific type is connected, you can use the [Mcs.Usb.CMcsUsbListNet](#) class to get a list of available devices:

```
CMcsUsbListNet usblast = new CMcsUsbListNet(DeviceEnumNet.MCS_DEVICE_USB);  
var entry = usblast.GetUsbListEntry((uint)0);  
CMcsUsbNet device = new CMcsUsbNet();  
device.Connect(entry);
```

After you are finished with the device, you can disconnect the device object from the device by:

```
device.Disconnect();
```

2 Device Classes

- For FluidControl device see [MCS FluidControl](#)
- For SW2TO64 device see [MCS-USB-Sw2to64](#)
- For TCx device see [Mcs.Usb.CTcxDeviceNet](#)

2.1 The MCS FluidControl Device

2.1.1 Introduction

The FluidControl Device can control up to 24 valves. The nominal voltage is 24V.

8 TTL level digital output ports are available and 8 TTL inputs can be read in.

The device has 8 ADC inputs with a range from 0V to 3.3V.

2.1.2 Access to the FluidControl device

For connecting to a FluidControl device see [Connecting to an MCS device.*](#)

```
CFluidControlDevice* m_dacq;  
m_fluidcontrol = new CFluidControlDevice;  
status = m_fluidcontrol->Connect();
```

The valves are controlled with the `CFluidControlDevice::SetValve` call. The argument given is a bit pattern of all valves which should be open.

The digital outputs can be controlled with the `CFluidControlDevice::SetDigout` call. Again, a bit pattern of all digital output pins which should be set to a logic high level is given as an argument.

The current state of the valves and the digital outputs can be read back with the `CFluidControlDevice::GetValve` and `CFluidControlDevice::GetDigout`

The command to read an ADC-Channel is `CFluidControlDevice::GetAdc`. Here the channelnummer which should be read in is given as an argument and the return value is the current Adc level.

The state of the digital inputs is read with the `CFluidControlDevice::GetDigin` call. Here the return value is the bit pattern of the digital inputs.

The connection to the device is closed with the `CFluidControlDevice::Disconnect` call.

2.2 MCS-USB-Sw2to64 device

The class [Mcs.Usb.CSw2to64DeviceNet](#) controls the setting of the switches in the MCS-USB-Sw2to64 device.

First construct an object of the class:

```
CSw2to64DeviceNet device = new CSw2to64DeviceNet();
```

For connecting to an MCS-USB-Sw2to64 device see [Connecting to an MCS device](#).

To get the number of channels the device handles:

```
int number = device.GetNumber();
```

Set all channel switches at once:

```
byte z = 1;
byte[] pattern = new byte[number];
for(int i = 0; i < number; i++)
{
    pattern[i] = z; // pattern you want to switch this channel to
}
device.SetChannels(pattern);
```

Get all channel switches at once:

```
byte[] pattern = device.GetChannels();
```

Set one channel switch:

```
ushort index = 10;
byte pattern = 1;
device.SetChannel(index, pattern)
```

Get one channel switch:

```
ushort index = 10;
byte pattern = device.GetChannel(index);
```

3 Function Classes

- [Mcs.Usb.CCMOSMea_FunctionNet](#)
- [Mcs.Usb.CDacCalibrationFunctionNet](#)
- [Mcs.Usb.CDigOutStimulatorFunctionNet](#)
- [Mcs.Usb.CGrapheneFunctionNet](#)
- [Mcs.Usb.CIntanMea_FunctionNet](#)
- [Mcs.Usb.CInterfaceboard2FunctionNet](#)
- [Mcs.Usb.CInterfaceboardFunctionNet](#)
- [Mcs.Usb.CMcsBus_MotorControlNet](#)
- [Mcs.Usb.CMcsBus_VoltageModeNet](#)
- [Mcs.Usb.CMcsBus_AxisParametersNet](#)
- [Mcs.Usb.CMcsBus_SensorNet](#)
- [Mcs.Usb.CMcsBus_TempSensorNet](#)
- [Mcs.Usb.CMcsBus_ExtensionNet](#)
- [Mcs.Usb.CMcsBus_FYIExtensionNet](#)
- [Mcs.Usb.CMcsUsbDeviceStatePushFunctionNet](#)
- [Mcs.Usb.CMEA2100x256FunctionNet](#)

- [Mcs.Usb.CMeaAudioFunctionNet](#)
- [Mcs.Usb.CMeaDigitalDataFunctionNet](#)
- [Mcs.Usb.CMeaFeedbackFunctionNet](#)
- [Mcs.Usb.CMeFunctionNet](#)
- [Mcs.Usb.CMultiwellCallbackFunctionNet](#)
- [Mcs.Usb.CMultiwellOptoStimFunctionNet](#)
- [Mcs.Usb.CPPCFunctionNet](#)
- [Mcs.Usb.CPPS_FunctionNet](#)
- [Mcs.Usb.CPulseGeneratorFunctionNet](#)
- [Mcs.Usb.CRFFunctionNet](#)
- [Mcs.Usb.CRobo_FYITemp_FunctionNet](#)
- [Mcs.Usb.CRobo_FYIProgram_FunctionNet](#)
- [Mcs.Usb.CSCUFunctionNet](#)
- [Mcs.Usb.CStimulusFunctionNet](#)
- [Mcs.Usb.CTEERFunctionNet](#)
- [Mcs.Usb.CW2100_FunctionNet](#)
- [Mcs.Usb.CW2100_StimulatorFunctionNet](#)
- [Mcs.Usb.CWarnerUssingFunctionNet](#)
- [Mcs.Usb.CWarnerValveControllerDeviceTesterFunctionNet](#)
- [Mcs.Usb.CWClassicFunctionNet](#)
- [Mcs.Usb.CWirelessBaseFunctionNet](#)

4 Data ACQuisition (DACQ) Devices

There are different device types of (MEA) data acquisition (DACQ) devices. All of them are supported by this class.

This library does **not** support the writing of the MCD (MC_Rack), MSRD (Multi Channel Experimenter) or HDF5 file format!

The class [Mcs.Usb.CMeaDeviceNet](#) is the base class for DACQ devices.

The base class [Mcs.Usb.CMeaDeviceNet](#) constructs actually the underlying classes for USB-MEA devices ([Mcs.Usb.CMeaUSBDeviceNet](#)).

```
CMeaDeviceNet device = new CMeaDeviceNet(McsBusTypeEnumNet.MCS_USB_BUS, OnChannelData, OnError);
```

For connecting to a DACQ device see [Connecting to an MCS device](#).

Get the number of available analog hardware channels and set the number of channels to the maximum.

```
int hwchannels;  
device.HWInfo().GetNumberOfHWADCCChannels(out hwchannels);  
device.SetNumberOfChannels(hwchannels);  
int samplingrate = 1000;  
device.SetSamplerate(samplingrate, 1, 0);  
device.EnableDigitalIn(true, 0);
```

Get the layout to know how the data look like that you receive

```
int ana, digi, che, tim, block;
device.GetChannelLayout(out ana, out digi, out che, out tim, out block);
```

For the [Mcs.Usb.OnChannelData](#) callback function you have to provide a definition of the channels you want to receive.

```
bool[] selChannels = new bool[block];
for (int i = 0; i < block; i++)
{
    selChannels[i] = true; // With true channel i is selected
    // selChannels[i] = false; // With false the channel i is deselected
}
channelblocksize = samplingrate / 10;
// queue size and threshold should be selected carefully
device.SetSelectedChannels(selChannels, 10 * channelblocksize, channelblocksize);
```

The [Mcs.Usb.OnChannelData](#) callback function gets a callback for each channelblock that is defined. In this example a callback for each channel.

```
void OnChannelData(CMcsUsbDacqNet d, int cbHandle, int numSamples)
{
    int size_ret;
    ushort[] channeldata = device.ChannelBlock_ReadFramesUI16(CbHandle, numSamples, out size_ret);
}
void OnError(String msg, int info)
{
    MessageBox.Show("Mea Device Error: " + msg);
}
```

see MEA_Recording in the Examples directory.

5 The MCS Robo Device

5.1 Introduction

Up to now two MCS devices exist that base on the Robo platform.

- The MCS Roboinject device is controlled by the [Mcs.Usb.CRoboInjectDeviceNet](#) class.
- The MCS Roboocyte2 device is controlled by the [Mcs.Usb.CRoboocyte2DeviceNet](#) class.

Both classes are derived from [Mcs.Usb.CRoboDeviceNet](#)

6 STG200x & STG400x STimulus Generator

6.1 Introduction

The STG200x & STG400x Series Stimulus Generators have two distinct modes of operation, the [Download mode](#) and the [Streaming mode](#).

6.2 Download mode

The Download mode is the "classic" mode of operation, as used by the MC Stimulus software. In this mode, one or multiple waveforms are defined in PC memory and downloaded to the STG. The waveforms are stored in STG device onboard memory and can be sent to the analog and sync outputs once or multiple times. The STG can operate independantly from the PC (without computer connection) after the download. Output is triggered either by the front panel start/stop button, the digital trigger inputs or under software control.

In the Download mode, there are up to eight independent triggers available (depending on the device). The user can assign each of the analog outputs and sync (digital) outputs to any of the triggers.

The analog output waveform is stored sample by sample in the STG memory. To reduce memory usage, this data can be compressed: whenever a given output value is to be held for more than one sample period, it has only to be given once. The user can define the number of sample periods for that a pattern should remain active. Compression is done for each channel independantly of the others, thus the algorithm to compress the data is very easy to implement.

A new feature of the Download mode is the segmentation of the STG memory. The onboard memory can be divided into up to 100 segments. Each segment can hold its own waveform pattern. Under software control, the user can switch between the defined segments within milliseconds. Another option is to use the four trigger inputs to select between four predefined segments. This option is accessible from the MC_Stimulus Software as the "Multi-File mode", and can start each of up to four defined waveforms within microseconds. This feature allows a predefined flexible response (feedback) to recorded data.

[Mcs.Usb.CStg200xDownloadNet](#) is the class for using the STG in download mode.

6.2.1 Memory Layout and Trigger Setup

The class to be used for the Download mode is [Mcs.Usb.CStg200xDownloadNet](#), which is derived from [Mcs.Usb.CStg200xBasicNet](#). You can add a poll handler delegate ([Mcs.Usb.OnStg200xPollStatus](#)) to the constructor [Mcs.Usb.CStg200xDownloadNet](#).

For connecting to an STG see [Connecting to an MCS device](#).

To use the Download mode, the memory layout of the STG200x can be set up, if the default is not sufficient. The total amount of memory available in the STG is obtained by the [Mcs.Usb.CStg200xDownloadNet.GetTotalMemory](#) call. With [Mcs.Usb.CStg200xDownloadNet.SendSegmentDefine](#) the segment sizes are assigned.

```
uint32_t memory = device.GetTotalMemory(); // obtain total memory available
uint[] segmentmemory = new uint[2];      // each segments has half of total memory
segmentmemory[0] = memory / 2;
segmentmemory[1] = memory / 2;
device.SendSegmentDefine(segmentmemory); // setup the STG
```

Next, for each segment, one has to assign the amount of memory to be used for each channel and sync output. This is done by [Mcs.Usb.CStg200xDownloadBasicNet.SetCapacity](#). Its arguments contain a list of memory sizes, with one entry per channel and one entry per sync output. Again, the total memory assigned to the channels and sync outputs must not exceed the memory assigned to the segment.

```
uint32_t nchannels = device.GetNumberOfAnalogChannels();
uint32_t nsync = device.GetNumberOfSyncoutChannels();
uint[] channel_cap = new uint[nchannels];
uint[] syncout_cap = new uint[nsync];
for (int i = 0; i < 2; i++) // for each segment
{
    device.SendSegmentSelect((uint32_t)i); // switch to segment
    uint32_t segment_mem = device.GetMemory(); // get memory available in this segment
    for(int j = 0; j < nchannels; j++)
    {
        channel_cap[j] = segment_mem/(nchannels+nsync); // devide memory amount to all channels
    }
    for(int j = 0; j < nsync; j++)
    {
        syncout_cap[j] = segment_mem/(nchannels+nsync); // and all sync outs.
    }
    device.SetCapacity(channel_cap, syncout_cap); // define memory for current segment
}
```

```
}
```

Before the STG can start, the trigger has to be configured. This is done by the [Mcs.Usb.CStg200xDownloadNet.SetupTrigger](#) call. Its arguments are a list of channelmaps, syncoutmaps and repeats, one for each of the four available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and syncout 1 to trigger 1 and channel 3 to trigger 2 use:

```
uint32_t TriggerInputs = device.GetNumberOfTriggerInputs();
uint[] channelmap = new uint[TriggerInputs];
uint[] syncoutmap = new uint[TriggerInputs];
uint[] repeat = new uint[TriggerInputs];
for (int i = 0; i < TriggerInputs; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    repeat[i] = 0;
}
// Trigger 0
channelmap[0] = 1; // Channel 1
syncoutmap[0] = 1; // Syncout 1
repeat[0] = 0; // forever
// Trigger 1
channelmap[1] = 4; // Channel 3
device.SetupTrigger(channelmap, syncoutmap, repeat);
```

For the STG400x series you have to set the output mode of the channels. [Mcs.Usb.CStg200xDownloadNet.SetVoltageMode](#) interprets the values as voltages. [Mcs.Usb.CStg200xDownloadNet.SetCurrentMode](#) as currents.

```
// Only meaningful for STG400x
device.SetVoltageMode();
```

For each segment, data can be sent to each of the defined channels and sync outputs using the [Mcs.Usb.CStg200xDownloadNet.SendChannelData](#) and [Mcs.Usb.CStg200xDownloadNet.SendSyncData](#) calls. channeldata and syncdata are a list of analog and digital samples as a list of two byte values (unsigned short). Multiple calls to [Mcs.Usb.CStg200xDownloadNet.SendChannelData](#) and [Mcs.Usb.CStg200xDownloadNet.SendSyncData](#) to the same channel append data to that channel.

If the Multi-File mode of the STG is enabled using the [Mcs.Usb.CStg200xDownloadNet.EnableMultiFileMode](#) call, the four trigger inputs are used to switch between four segments. A hardware trigger signal (TTL) on trigger input 1 selects the first segment and starts all pulses in this segment. Thus with the Multi-File mode, one can predefine four stimulus patterns and switch between them without a connection to the PC.

The STG200x series has an analog resolution of 13 bits, thus the analog data contains the information in bits 0 to 12 of each sample. Bits 13 to 15 have to be 0.

```
int DACResolution = device.GetDACResolution();
// Data for Channel 0
{
    device.ClearChannelData(0);
    double factor = 0.1;
    const int l = 1000;
    ushort[] pData = new ushort[l];
    UInt64_t[] tData = new UInt64_t[l];
    for (int i = 0; i < l; i++)
    {
        // calculate Sin-Wave
        double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
            Math.Sin(2.0 * (double)i * Math.PI / (double)l);
        // calculate sign
        pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
            (int)Math.Pow(2, DACResolution - 1));
        tData[i] = (UInt64_t)20; // duration in µs
    }
    device.SendChannelData(0, pData, tData);
}
// Data for Channel 3
{
    device.ClearChannelData(2);
    double factor = 0.1;
    const int l = 700;
    // without compression
    ushort[] pData = new ushort[l];
    UInt64_t[] tData = new UInt64_t[l];
    for (int i = 0; i < l; i++)
    {
        // calculate Sin-Wave
        double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
            Math.Sin(2.0 * (double)i * Math.PI / (double)l);
        // calculate sign
```

```

        pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
            (int)Math.Pow(2, DACResolution - 1));
        tData[i] = (UInt64_t)20; // duration in µs
    }
    device.SendChannelData(2, pData, tData);
}
// Data for Sync 0
{
    device.ClearSyncData(0);
    ushort[] pData = new ushort[1000];
    UInt64_t[] tData = new UInt64_t[1000];
    for (int i = 0; i < 1000; i++)
    {
        pData[i] = (ushort)(i&1);
        tData[i] = 20;
    }
    device.SendSyncData(0, pData, tData);
}

```

Start the trigger by pushing the front button or by software

```

// Start Trigger 1 and 2
device.SendStart(1 + 2); // Trigger 1 und 2

```

see the StgDownloadExampleNet in the example directory.

6.3 Streaming mode

The other mode of operation is the Streaming mode. Here the analog output is sent to the STG device in "real time". The PC has to be connected to the STG all the time. The data that is sent to the analog output is downloaded from the PC to the STG on the fly.

The Streaming mode is useful for applications where flexible feedback is needed as well for applications where very long waveforms which are not repeated (such as white noise) are used.

The Streaming mode works by use of two ring buffers which hold data. One is in PC memory and managed by the DLL, and one is in on-board STG memory. Data is transferred from PC memory to the STG via the USB bus in time slices of one millisecond.

The user can define both the size of the ring buffer in DLL memory and in the STG memory. Once the Streaming mode is started, the STG request data from the PC. The data rate from PC to STG is variable and controlled by the STG. The STG request data from the PC at a rate to keep its internal ringbuffer at about half full.

It is the responsibility of the user to keep the ring buffer in the memory of the PC filled, so the DLL can supply sufficient data to the STG. To do so, the Windows DLL allows to define a "callback" function which is called whenever new data is needed, or more precise, as soon as the ring buffer in the memory of the PC falls below the user defined threshold.

Small buffers have the advantage of a low latency between data generation in the callback function and its output as a analog signal from the STG. However for low latency to work, the user-written callback function has to be fast and to produce a steady flow of data.

In the Streaming mode, all triggers are available as well. Each of the eight analog and sync outputs can be assigned to one of the triggers.

The output rate is user defined with a maximum of 50 kHz

Mcs.Usb.CStg200xStreamingNet is the class for using the STG in streaming mode.

6.3.1 Memory Layout and Trigger Setup

With the constructor for `Mcs.Usb.CStg200xStreamingNet.CStg200xStreamingNet`, the name of the callback function for the data handler is provided. The data handler function is called automatically, whenever the STG needs new data. This data is first written to a ring buffer in the memory of the PC. The size for this ring buffer is defined as first argument in the constructor. The user provided delegate gets the trigger number which needs new data as argument

```
CStg200xStreamingNet device = new CStg200xStreamingNet(10000, dataHandler, errorHandler);
```

The callback function, which is defined in the constructor, is called whenever the STG needs new data for a trigger, or more precise, whenever the ring buffer in PC memory falls below the defined threshold.

The user can query the amount of space available for queuing by use of the `Mcs.Usb.CStg200xStreamingNet.GetDataQueueSpace` call. Its return value is the number of samples that can be send to the STG.

User code is required to fill an array analog and sync out data, sample by sample for up to the maximum number of samples as obtained by `Mcs.Usb.CStg200xStreamingNet.GetDataQueueSpace` or `Mcs.Usb.CStg200xStreamingNet.GetSyncoutQueueSpace`.

The values for the analog outputs are 16 bits signed integers. The lower bits are truncated according to the resolution of the STG. This behaviour is different to the behaviour in [download mode](#).

Note: Compression as described in the [download mode](#) can NOT be used for the streaming mode.

The new data is sent to the STG by using the `Mcs.Usb.CStg200xStreamingNet.EnqueueData` call.

```
void dataHandler(uint32_t trigger)
{
    double factor = 1;
    if (trigger == 0) // Callback for Trigger 1
    {
        // Handle Channel 1
        uint32_t channel = 0;
        for ( ; ; )
        {
            uint32_t space = device.GetDataQueueSpace(channel);
            if (space < 1000)
                break;
            short[] data = new short[1000];
            for (int i = 0; i < 1000; i++)
            {
                // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
                double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) *
                    Math.Sin(2.0 * (double)i * Math.PI / (double)1000);
                data[i] = (short)sin;
            }
            uint32_t enqueued = device.EnqueueData(channel, data);
        }
    }
    // Handle Channel 3
    uint32_t channel = 2;
    for ( ; ; )
    {
        uint32_t space = device.GetDataQueueSpace(channel);
        if (space < 700)
            break;
        short[] data = new short[700];
        for (int i = 0; i < 700; i++)
        {
            // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
            double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) *
                Math.Sin(2.0 * (double)i * Math.PI / (double)700);
            data[i] = (short)sin;
        }
        uint32_t enqueued = device.EnqueueData(channel, data);
    }
}
// Handle Syncout 1
uint32_t channel = 0;
for ( ; ; )
{
    uint32_t space = device.GetSyncoutQueueSpace(channel);
    if (space < 1000)
        break;
    ushort[] data = new ushort[1000];
    for (int i = 0; i < 1000; i++)
```

```

        data[i] = (ushort)(i & 1);
        uint32_t enqueued = device.EnqueueSyncout(channel, data);
    }
}
}
}
void errorHandler()
{
}

```

For connecting to an STG device see [Connecting to an MCS device](#).

With enabling or disabling the continuous mode it can be selected how the STG handles an "out of data" situation.

When `Mcs.Usb.CStg200xStreamingNet.EnableContinuousMode` is used, the STG does not stop when it runs out of data, but it keeps running and sends a zero voltage to its outputs.

When `Mcs.Usb.CStg200xStreamingNet.DisableContinuousMode` is used, the STG stops when it runs out of data. It has to be retriggered to resume the output.

```
device.EnableContinuousMode();
```

`Mcs.Usb.CStg200xStreamingNet.SetOutputRate` is used to set the sampling rate.

```
device.SetOutputRate(50000);
```

To use the Streaming mode, the memory layout of the STG has to be set up. To total amount of memory available in the STG is obtained by the `Mcs.Usb.CStg200xStreamingNet.GetTotalMemory` call.

This memory can be assigned to four ring buffers (one per trigger) which buffer the data received from the PC via USB cable. This is done with the `CStg200xStreaming::SetCapacity` call. The total amount of memory must not exceed the total memory size as obtained by `Mcs.Usb.CStg200xStreamingNet.GetTotalMemory`.

This internal ring buffer is crucial for proper operation of the Streaming mode. The size of the ring buffer determines the latency of the Streaming mode. The firmware of the STG requests data from the PC in order to keep the ring buffer about half full. Thus the average latency is:

$$\text{latency} = (\text{ringbuffersize in bytes}/4) / \text{output rate}$$

If the ring buffer size is too big, the latency of the STG might be too long. If the ring buffer size is too low, an overflow or underflow of data in the STG ringbuffer might occur, resulting in data jumps of the output signals or the "out of data" situation described earlier.

The following example divides the total memory equally among the four triggers:

```

uint32_t dwMemory = device.GetTotalMemory(); // obtain total memory available
uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] stg_triggercapacity = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
    stg_triggercapacity[i] = dwMemory / ntrigger;
device.SetCapacity(stg_triggercapacity); // setup the STG

```

or fixed memory sizes:

```

uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] stg_triggercapacity = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
    stg_triggercapacity[i] = 50000;
device.SetCapacity(stg_triggercapacity);

```

Before the STG can start, the trigger has to be configured. This is done by the `Mcs.Usb.CStg200xStreamingNet.SetupTrigger` call. Its arguments are a list of channelmaps, syncoutmaps, digoutmap, autostart and callback←_threshold, with one entry for each of the available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and 3 and syncout 1 to trigger 1 use:

```

uint32_t ntrigger = device.GetNumberOfTriggerInputs(); // obtain number of triggers in this STG
uint[] channelmap = new uint[ntrigger];
uint[] syncoutmap = new uint[ntrigger];
uint[] digoutmap = new uint[ntrigger];
uint[] autostart = new uint[ntrigger];
uint[] callback_threshold = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
{
    channelmap[i] = 0;

```



```

    syncoutmap[i] = 0;
    digoutmap[i] = 0;
    autostart[i] = 0;
    callback_threshold[i] = 0;
}
channelmap[0] = 0x1 + 0x4; // Channel 1 und Channel 3 to Trigger 1
syncoutmap[0] = 0x1; // Syncout 1 to Trigger 1
autostart[0] = 1;
callback_threshold[0] = 50; // 50% of buffer size
device.SetupTrigger(channelmap, syncoutmap, digoutmap, autostart, callback_threshold);
device.StartLoop();
System.Threading.Thread.Sleep(1000); // Give StartLoop some time

```

Start Trigger by pushing the front button or by Software

```
device.SendStart(1);
```

see the StgStreamingExampleNet in the example directory.

7 Namespace Index

7.1 Namespace List

Here is a list of all namespaces with brief descriptions:

Mcs	22
Mcs::Usb	22

8 Hierarchical Index

8.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

CW2100_FunctionNet::AudioChannelsNet	92
BatteryState	92
CCreateFilterNet	110
BesselFilterHighPassNet	93
BesselFilterLowPassNet	93
ButterworthFilterHighPassNet	94
ButterworthFilterLowPassNet	95
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet >	120
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumNet >	120
CDeviceGroupChannelInfoNet	119
CDeviceGroupChannelInfoTemplateNet< int >	120
CDeviceGroupChannelInfoGenericNet	118

CDeviceGroupChannelInfoTemplateNet< MEA2100_256DacqGroupChannelEnumNet >	120
CDeviceGroupChannelInfoMEA2100_256Net	119
CDeviceGroupChannelInfoTemplateNet< SCUDacqGroupChannelEnumNet >	120
CDeviceGroupChannelInfoSCUNet	120
CDeviceGroupChannelInfoTemplateNet< W2100DacqGroupChannelEnumNet >	120
CDeviceGroupChannelInfoW2100Net	121
CFilterCoefficientsNet	129
CFilterPropertyNet	135
CMcsUsbDacqNet::CHWInfo	179
CMcsUsbFunctionNet	298
CDacqGroupChannelSelectionTemplateNet< W2100DacqGroupChannelEnumNet, W2100DacqGroupChannelEnum, CDeviceGroupChannelInfoW2100Net >	115
CW2100DacqGroupChannelSelectionNet	611
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumNet, DacqGroupChannelEnum, CDeviceGroupChannelInfoNet >	115
CDacqGroupChannelSelectionNet	115
CDacqGroupChannelSelectionTemplateNet< int, int, CDeviceGroupChannelInfoGenericNet >	115
CDacqGroupChannelGenericSelectionNet	114
CDacqGroupChannelSelectionTemplateNet< SCUDacqGroupChannelEnumNet, SCUDacqGroupChannelEnum, CDeviceGroupChannelInfoSCUNet >	115
CSCUDacqGroupChannelSelectionNet	491
CDacqGroupChannelSelectionTemplateNet< MEA2100_256DacqGroupChannelEnumNet, MEA2100_256DacqGroupChannelEnum, CDeviceGroupChannelInfoMEA2100_256Net >	115
CMEA2100_256DacqGroupChannelSelectionNet	329
CCMOSMea_FunctionNet	96
CDacCalibrationFunctionNet	112
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >	115
CDigOutStimulatorFunctionNet	122
CFilterConfigurationNet	131
CFilterConfigurationRegisterNet	133
CGrapheneFunctionNet	166
CIntanMea_FunctionNet	182
CInterfaceboardFunctionNet	186

CInterfaceboard2FunctionNet	184
CMEA2100x256FunctionNet	329
CMcsBusNet	236
CMcsBus_AxisParametersNet	197
CMcsBus_ExtensionNet	199
CMcsBus_FYIExtensionNet	200
CMcsBus_MotorControlNet	202
CMcsBus_SensorNet	219
CMcsBus_TempSensorNet	229
CMcsBus_VoltageModeNet	231
CMcsUsbDeviceStatePushFunctionNet	287
CMultiwellCallbackFunctionNet	372
CSCUFunctionNet	492
CMeFunctionNet	363
CMeaAudioFunctionNet	331
CMeaDigitalDataFunctionNet	351
CMeaFeedbackFunctionNet	353
CMultiwellOptoStimFunctionNet	382
CPPCFunctionNet	414
CPPS_FunctionNet	424
CProgramPressureCurveNet	429
CPulseGeneratorFunctionNet	430
CRFFunctionNet	438
CRobo_FYIProgram_FunctionNet	443
CRobo_FYITemp_FunctionNet	444
CStimulusFunctionNet	556
CTEERFunctionNet	582
CUsbDeviceConfigurationFunctionNet	594
CW2100_StimulatorFunctionNet	605
CWarnerUssingFunctionNet	613
CWarnerValveControllerDeviceTesterFunctionNet	654
CWirelessBaseFunctionNet	662

CW2100_FunctionNet	597
CWClassicFunctionNet	657
CMcsUsbFunctionPointerContainer	299
CMcsUsbListEntryNet	299
CMcsUsbListNet	305
CMcsUsbNet	308
CExternDTesterDeviceNet	127
CFluidControlDeviceNet	136
CGenericDevelopDeviceNet	145
CGilsonDeviceNet	163
CMcsUsbDacqNet	240
CMeaDeviceNet	343
CMeaUSBDeviceNet	361
CCMOSMeaDeviceNet	107
CGrapheneASICDeviceNet	165
CHLADacqNet	178
CLIH3DeviceNet	188
CMultiwellDeviceNet	374
CWarnerUssingDeviceNet	612
COctoPotDeviceNet	387
CRoboDacqNet	446
CMcsUsbDeviceStatePushNet	288
CWarnerValveControllerDeviceNet	629
CMcsUsbFactoryNet	289
CMeaCleanDeviceNet	334
CMeaCoatDeviceNet	338
CMeaImpedanceDeviceNet	357
CMeaSwitchDeviceNet	360
CChannelTestDeviceNet	95
CMultiBatteryChargerDeviceNet	364
CNF_GenDeviceNet	386
COkuvisionStimulatorDeviceNet	391

CPPCDeviceNet	413
CPPS_DeviceNet	423
CPathIdentDeviceNet	396
CPedoterDeviceNet	397
CPeristalticPumpDeviceNet	398
CPgaDeviceNet	399
CPositionIIDeviceNet	401
CPositionImpDeviceNet	410
CRadioControlledDevicesNet	433
CRetinaLedDeviceNet	436
CRoboDeviceNet	461
CEncapsulatorDeviceNet	126
CFYIDeviceNet	143
CHLADeviceNet	178
CHiClampDeviceNet	177
CMeasureTableDeviceNet	359
CPPSDeviceNet	428
CPatchServerDeviceNet	395
CRoboStatorDeviceNet	482
CRoboInjectDeviceNet	480
CRoboocyte2DeviceNet	481
CTEERMachineDeviceNet	593
CRoboFluidDeviceNet	477
CSafeISDeviceNet	488
CSerialPortNet	507
CStg200xBasicNet	508
CStg200xDownloadBasicNet	542
CStg200xDownloadNet	549
CSw2to64DeviceNet	567
CTcxDeviceNet	569
CMcsUsbPointerContainer	329
CCMOSMeaDeviceNet::CRegionOfInterestRect	435

CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	596
DeviceIdNet	663
DigitalSource< digitalsourceenum >	665
DigitalSourceGeneral	666
DriverVersionNet	667
Exception	
CUsbExceptionNet	595
FirmwareDestinationNames	673
HeadstageIDTypeObject	679
HeadStageIDTypeState	681
IComparable	
HeadStageIDType	676
mkfilterNet	682
CRoboDeviceNet::RoboMainLowLevelCommands	685
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	692
CFilterCoefficientsNet::s_FilterAttributesNet	692
CMeaAudioFunctionNet::s_setaudionet	693
CStimulusFunctionNet::SidebandData	694
StgStatusNet	695
stgstreaming	
CStg200xBasicNet	508
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	695
usbSetupPacket_t	697
W2100_StimulusParametersNet	697

9 Class Index

9.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

CW2100_FunctionNet::AudioChannelsNet	92
BatteryState	92
BesselFilterHighPassNet	93
BesselFilterLowPassNet	93

ButterworthFilterHighPassNet	94
ButterworthFilterLowPassNet	95
CChannelTestDeviceNet	95
CCMOSMea_FunctionNet	96
CCMOSMeaDeviceNet	107
CCreateFilterNet	110
CDacCalibrationFunctionNet 112	
CDacqGroupChannelGenericSelectionNet	114
CDacqGroupChannelSelectionNet	115
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplateNet > 115	
CDeviceGroupChannelInfoGenericNet	118
CDeviceGroupChannelInfoMEA2100_256Net	119
CDeviceGroupChannelInfoNet	119
CDeviceGroupChannelInfoSCUNet	120
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet >	120
CDeviceGroupChannelInfoW2100Net	121
CDigOutStimulatorFunctionNet CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class	122
CEncapsulatorDeviceNet CEncapsulatorDeviceNet is the to control the MCS HiClamp device	126
CExternDTesterDeviceNet CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)	127
CFilterCoefficientsNet	129
CFilterConfigurationNet	131
CFilterConfigurationRegisterNet	133
CFilterPropertyNet	135
CFluidControlDeviceNet CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device	136
CFYIDeviceNet CFYIDeviceNet is the class to control the MCS FYI device	143
CGenericDevelopDeviceNet CGenericDevelopDeviceNet is the class to use during development of a new device	145
CGilsonDeviceNet CGilsonDeviceNet is the class to control a Gilson device	163

CGrapheneASICDeviceNet	165
CGrapheneFunctionNet	
CGrapheneFunctionNet is the class to control Graphene device functions	166
CHiClampDeviceNet	
CHiClampDeviceNet is the to control the MCS HiClamp device	177
CHLADacqNet	178
CHLADeviceNet	
CHLADeviceNet is the to control the MCS HLA device	178
CMcsUsbDacqNet::CHWInfo	
Class to provide hardware information about the device	179
CIntanMea_FunctionNet	182
CInterfaceboard2FunctionNet	
CInterfaceboard2FunctionNet is the class to control the Interfaceboard	184
CInterfaceboardFunctionNet	
CInterfaceboardFunctionNet is the class to control the Interfaceboard	186
CLIH3DeviceNet	
CLIH3DeviceNet is the class to access the HEKA LIH3 device	188
CMcsBus_AxisParametersNet	197
CMcsBus_ExtensionNet	199
CMcsBus_FYIExtensionNet	200
CMcsBus_MotorControlNet	202
CMcsBus_SensorNet	219
CMcsBus_TempSensorNet	229
CMcsBus_VoltageModeNet	231
CMcsBusNet	236
CMcsUsbDacqNet	
Base class for data acquisition devices	240
CMcsUsbDeviceStatePushFunctionNet	287
CMcsUsbDeviceStatePushNet	288
CMcsUsbFactoryNet	289
CMcsUsbFunctionNet	298
CMcsUsbFunctionPointerContainer	299
CMcsUsbListEntryNet	
McsUsbListEntryNet identifies a connected device	299
CMcsUsbListNet	
Class to handle a list of connected MCS USB devices	305

CMcsUsbNet	
Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class	308
CMcsUsbPointerContainer	329
CMEA2100_256DacqGroupChannelSelectionNet	329
CMEA2100x256FunctionNet	
CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "↔ Stg200xNet.h" to resolve documentation reference	329
CMeaAudioFunctionNet	331
CMeaCleanDeviceNet	
CMeaCleanDeviceNet is the class to access the MEA Clean device	334
CMeaCoatDeviceNet	
CMeaCoatDeviceNet is the class to access the MEA Coat device	338
CMeaDeviceNet	
Base class for MEA data acquisition devices	343
CMeaDigitalDataFunctionNet	351
CMeaFeedbackFunctionNet	353
CMeaImpedanceDeviceNet	357
CMeasureTableDeviceNet	
CMeasureTableDeviceNet is the to control the MCS HLA device	359
CMeaSwitchDeviceNet	
The class to control the USB-MEA-Switch	360
CMeaUSBDeviceNet	
Class for data acquisition via ME and MEA USB amplifiers	361
CMeFunctionNet	363
CMultiBatteryChargerDeviceNet	
CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device	364
CMultiwellCallbackFunctionNet	
CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator	372
CMultiwellDeviceNet	
CMultiwellDeviceNet is the class to access the Multiwell device	374
CMultiwellOptoStimFunctionNet	
CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device	382
CNF_GenDeviceNet	386
COctoPotDeviceNet	387
COKuvisionStimulatorDeviceNet	391
CPatchServerDeviceNet	
CPatchServerDeviceNet is the class to control the MCS PatchServer device	395

CPathIdentDeviceNet	396
CPedoterDeviceNet	397
CPeristalticPumpDeviceNet	
CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump	398
CPgaDeviceNet	399
CPositionIIDeviceNet	
CPositionIIDeviceNet is the class to control PositionII devices	401
CPositionImpDeviceNet	
CPositionImpDeviceNet is the class to access the Position/Imp devices	410
CPPCDeviceNet	413
CPPCFunctionNet	
CPPCFunctionNet is the class to access the PPC (high precision Patch Peristaltic patch Pump	414
CPPS_DeviceNet	423
CPPS_FunctionNet	424
CPPSDeviceNet	
CPPS4plus1DeviceNet is the to control the MCS HLA device	428
CProgramPressureCurveNet	
CProgramPressureCurveNet is the class to program pressure curves	429
CPulseGeneratorFunctionNet	
CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking	430
CRadioControlledDevicesNet	433
CCMOSMeaDeviceNet::CRegionOfInterestRect	435
CRetinaLedDeviceNet	436
CRFFFunctionNet	
CRFFFunctionNet is the class to control RF devices	438
CRobo_FYIProgram_FunctionNet	443
CRobo_FYITemp_FunctionNet	444
CRoboDacqNet	446
CRoboDeviceNet	
CRoboDeviceNet is the base class for all Robo platform based devices	461
CRoboFluidDeviceNet	477
CRoboInjectDeviceNet	
CRoboInjectDeviceNet is the to control the MCS RoboInject device	480
CRoboocyte2DeviceNet	
CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device	481
CRoboStatorDeviceNet	482

CSafeISDeviceNet 488	
CSCUDacqGroupChannelSelectionNet	491
CSCUFunctionNet CSCUFunctionNet is the class to control the SCU device	492
CSerialPortNet	507
CStg200xBasicNet Base class for the Stg200x	508
CStg200xDownloadBasicNet CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device	542
CStg200xDownloadNet Main class for the STG download mode This class implements the STG download mode interface.	549
CStimulusFunctionNet	556
CSw2to64DeviceNet The class to control the MCS-USB-Sw2to64 device	567
CTcxDeviceNet Class to control a Temperature Controller (TCX)	569
CTEERFunctionNet CTEERFunctionNet is the class to control the TEER device	582
CTEERMachineDeviceNet	593
CUsbDeviceConfigurationFunctionNet CUsbDeviceConfigurationFunctionNet is the class to configure the USB firmware	594
CUsbExceptionNet Exception class that is thrown in case of an USB error	595
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	596
CW2100_FunctionNet	597
CW2100_StimulatorFunctionNet	605
CW2100DacqGroupChannelSelectionNet	611
CWarnerUssingDeviceNet CWarnerUssingDeviceNet is the class to control the Ussing device	612
CWarnerUssingFunctionNet CWarnerUssingFunctionNet is the class to control the Ussing device	613
CWarnerValveControllerDeviceNet CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller	629
CWarnerValveControllerDeviceTesterFunctionNet CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester	654
CWClassicFunctionNet	657

CWirelessBaseFunctionNet	662
DeviceIdNet	
Device Id	663
DigitalSource< digitalsourceenum >	665
DigitalSourceGeneral	666
DriverVersionNet	
Class gives firmware versions of the device's firmware destinations	667
FirmwareDestinationNames	673
HeadStageIDType	676
HeadstageIDTypeObject	679
HeadStageIDTypeState	681
mkfilterNet	682
CRoboDeviceNet::RoboMainLowLevelCommands	685
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	692
CFilterCoefficientsNet::s_FilterAttributesNet	692
CMeaAudioFunctionNet::s_setaudionet	693
CStimulusFunctionNet::SidebandData	694
StgStatusNet	695
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	695
usbSetupPacket_t	697
W2100_StimulusParametersNet	697

10 Namespace Documentation

10.1 Mcs Namespace Reference

Namespaces

- [Usb](#)

10.2 Mcs::Usb Namespace Reference

Classes

- class [CChannelTestDeviceNet](#)
- class [CCMOSMea_FunctionNet](#)
- class [CCMOSMeaDeviceNet](#)
- class [CCreateFilterNet](#)

- class [ButterworthFilterLowPassNet](#)
- class [ButterworthFilterHighPassNet](#)
- class [BesselFilterLowPassNet](#)
- class [BesselFilterHighPassNet](#)
- class [CDeviceGroupChannelInfoTemplateNet](#)
- class [CDeviceGroupChannelInfoGenericNet](#)
- class [CDeviceGroupChannelInfoNet](#)
- class [CDeviceGroupChannelInfoW2100Net](#)
- class [CDeviceGroupChannelInfoSCUNet](#)
- class [CDeviceGroupChannelInfoMEA2100_256Net](#)
- class [CDacqGroupChannelSelectionTemplateNet](#)
- class [CDacqGroupChannelGenericSelectionNet](#)
- class [CDacqGroupChannelSelectionNet](#)
- class [CW2100DacqGroupChannelSelectionNet](#)
- class [CSCUDacqGroupChannelSelectionNet](#)
- class [CMEA2100_256DacqGroupChannelSelectionNet](#)
- class [CDacCalibrationFunctionNet](#)
- class [CDigOutStimulatorFunctionNet](#)
 - [CDigOutStimulatorFunctionNet](#) is the class of the DigOut stimulator function class.*
- class [CExternDTesterDeviceNet](#)
 - [CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)*
- class [CGrapheneFunctionNet](#)
 - [CGrapheneFunctionNet](#) is the class to control Graphene device functions*
- class [CInterfaceboard2FunctionNet](#)
 - [CInterfaceboard2FunctionNet](#) is the class to control the Interfaceboard*
- class [CInterfaceboardFunctionNet](#)
 - [CInterfaceboardFunctionNet](#) is the class to control the Interfaceboard*
- class [CLIH3DeviceNet](#)
 - [CLIH3DeviceNet](#) is the class to access the HEKA LIH3 device.*
- class [CMEA2100x256FunctionNet](#)
 - [CMEA2100x256FunctionNet](#) is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference*
- class [CMeaCleanDeviceNet](#)
 - [CMeaCleanDeviceNet](#) is the class to access the MEA Clean device.*
- class [CMeaCoatDeviceNet](#)
 - [CMeaCoatDeviceNet](#) is the class to access the MEA Coat device.*
- class [CMultiBatteryChargerDeviceNet](#)
 - [CMultiBatteryChargerDeviceNet](#) is the class to access the MBC-08 device.*
- class [CMultiwellCallbackFunctionNet](#)
 - [CMultiwellCallbackFunctionNet](#) is the class to access the Multiwell-Mini-Stimulator*
- class [CMultiwellDeviceNet](#)
 - [CMultiwellDeviceNet](#) is the class to access the Multiwell device.*
- class [CMultiwellOptoStimFunctionNet](#)
 - [CMultiwellOptoStimFunctionNet](#) is the class to access the optical properties of the Multiwell Optostim device*
- class [CPedoterDeviceNet](#)
- class [CPositionIIDeviceNet](#)
 - [CPositionIIDeviceNet](#) is the class to control PositionII devices*
- class [CPositionImpDeviceNet](#)
 - [CPositionImpDeviceNet](#) is the class to access the Position/Imp devices*
- class [CPPCFunctionNet](#)
 - [CPPCFunctionNet](#) is the class to access the PPC (high precision Patch Peristaltic patch Pump*
- class [CPulseGeneratorFunctionNet](#)

- CPulseGeneratorFunctionNet* is the class to control the pulse generator for video tracking
- class [CRFFunctionNet](#)
 - CRFFunctionNet* is the class to control RF devices
- class [CSCUFunctionNet](#)
 - CSCUFunctionNet* is the class to control the SCU device
- class [CTEERFunctionNet](#)
 - CTEERFunctionNet* is the class to control the TEER device
- class [CUsbDeviceConfigurationFunctionNet](#)
 - CUsbDeviceConfigurationFunctionNet* is the class to configure the USB firmware
- class [CWarnerUssingDeviceNet](#)
 - CWarnerUssingDeviceNet* is the class to control the Ussing device
- class [CWarnerUssingFunctionNet](#)
 - CWarnerUssingFunctionNet* is the class to control the Ussing device
- class [CWarnerValveControllerDeviceNet](#)
 - CWarnerValveControllerDeviceNet* is the class to access the Warner Valve Controller
- class [CWarnerValveControllerDeviceTesterFunctionNet](#)
 - CWarnerValveControllerDeviceTesterFunctionNet* is the class to access the functions for the Warner Valve Controller Device Tester
- struct [DeviceIdNet](#)
 - Device Id.*
- class [CFilterCoefficientsNet](#)
- class [CFilterConfigurationNet](#)
- class [CFilterConfigurationRegisterNet](#)
- class [CFilterPropertyNet](#)
- class [CFluidControlDeviceNet](#)
 - CFluidControlDeviceNet* is the class to control MCS FluidControl (FCB and FCX) device.
- class [CGenericDevelopDeviceNet](#)
 - CGenericDevelopDeviceNet* is the class to use during development of a new device.
- class [CGilsonDeviceNet](#)
 - CGilsonDeviceNet* is the class to control a Gilson device.
- class [CGrapheneASICDeviceNet](#)
- class [CIntanMea_FunctionNet](#)
- class [CMcsBusNet](#)
- class [CMcsBus_MotorControlNet](#)
- class [CMcsBus_VoltageModeNet](#)
- class [CMcsBus_AxisParametersNet](#)
- class [CMcsBus_SensorNet](#)
- class [CMcsBus_TempSensorNet](#)
- class [CMcsBus_ExtensionNet](#)
- class [CMcsBus_FYIExtensionNet](#)
- class [CSerialPortNet](#)
- class [usbSetupPacket_t](#)
- class [CMcsUsbDeviceStatePushFunctionNet](#)
- class [CMcsUsbDeviceStatePushNet](#)
- class [CMcsUsbFactoryNet](#)
- class [CMcsUsbFunctionPointerContainer](#)
- class [CMcsUsbFunctionNet](#)
- class [CMcsUsbListEntryNet](#)
 - McsUsbListEntryNet* identifies a connected device.
- class [CMcsUsbListNet](#)
 - Class to handle a list of connected MCS USB devices.*
- class [CUsbExceptionNet](#)

Exception class that is thrown in case of an USB error.

- class [FirmwareDestinationNames](#)
- class [DriverVersionNet](#)

Class gives firmware versions of the device's firmware destinations.

- class [CMcsUsbPointerContainer](#)
- class [CMcsUsbNet](#)

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

- class [DigitalSourceGeneral](#)
- class [DigitalSource](#)
- class [StgStatusNet](#)
- class [CMeaAudioFunctionNet](#)
- class [CMeaDeviceNet](#)

Base class for MEA data acquisition devices.

- class [CMeaUSBDeviceNet](#)

Class for data acquisition via ME and MEA USB amplifiers

- class [CMeaDigitalDataFunctionNet](#)
- class [CMeaFeedbackFunctionNet](#)
- class [CMeaImpedanceDeviceNet](#)
- class [CMeaSwitchDeviceNet](#)

The class to control the USB-MEA-Switch.

- class [CMeFunctionNet](#)
- class [mkfilterNet](#)
- class [CNF_GenDeviceNet](#)
- class [COctoPotDeviceNet](#)
- class [COkuvisionStimulatorDeviceNet](#)
- class [CPathIdentDeviceNet](#)
- class [CPeristalticPumpDeviceNet](#)

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

- class [CPgaDeviceNet](#)
- class [CPPCDeviceNet](#)
- class [CPPS_DeviceNet](#)
- class [CPPS_FunctionNet](#)
- class [CProgramPressureCurveNet](#)

CProgramPressureCurveNet is the class to program pressure curves

- class [CRadioControlledDevicesNet](#)
- class [CRetinaLedDeviceNet](#)
- class [CRobo_FYITemp_FunctionNet](#)
- class [CRobo_FYIProgram_FunctionNet](#)
- class [CRoboDacqNet](#)
- class [CHLADacqNet](#)
- class [CRoboDeviceNet](#)

CRoboDeviceNet is the base class for all Robo platform based devices

- class [CRoboStatorDeviceNet](#)
- class [CRoboocyte2DeviceNet](#)

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

- class [CRoboInjectDeviceNet](#)

CRoboInjectDeviceNet is the to control the MCS RoboInject device

- class [CHiClampDeviceNet](#)

CHiClampDeviceNet is the to control the MCS HiClamp device

- class [CEncapsulatorDeviceNet](#)

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

- class [CHLADeviceNet](#)

- CHLADeviceNet* is the to control the MCS HLA device
- class [CPPSDeviceNet](#)
 - CPPS4plus1DeviceNet* is the to control the MCS HLA device
- class [CMeasureTableDeviceNet](#)
 - CMeasureTableDeviceNet* is the to control the MCS HLA device
- class [CFYIDeviceNet](#)
 - CFYIDeviceNet* is the class to control the MCS FYI device
- class [CPatchServerDeviceNet](#)
 - CPatchServerDeviceNet* is the class to control the MCS PatchServer device
- class [CTEERMachinDeviceNet](#)
- class [CRoboFluidDeviceNet](#)
- class [CSafeISDeviceNet](#)
- class [CStg200xDownloadNet](#)
 - Main class for the STG download mode This class implements the STG download mode interface.*
- class [CStg200xBasicNet](#)
 - Base class for the Stg200x.*
- class [CStg200xDownloadBasicNet](#)
 - CStg200xDownloadBasicNet* is the base class to control the download mode of the MCS STG device.
- class [CStimulusFunctionNet](#)
- class [CSw2to64DeviceNet](#)
 - The class to control the MCS-USB-Sw2to64 device.*
- class [CTcxDeviceNet](#)
 - Class to control a Temperature Controller (TCX)*
- class [CMcsUsbDacqNet](#)
 - Base class for data acquisition devices.*
- struct [W2100_StimulusParametersNet](#)
- class [HeadStageIDType](#)
- class [HeadStageIDTypeState](#)
- class [HeadstageIDTypeObject](#)
- class [BatteryState](#)
- class [CW2100_StimulatorFunctionNet](#)
- class [CW2100_FunctionNet](#)
- class [CWClassicFunctionNet](#)
- class [CWirelessBaseFunctionNet](#)

Enumerations

- enum class [enCMosMeaChipType](#) {
 - [unknown](#) = 0 ,
 - [nMos16LV](#) = 1 ,
 - [nMos32LV](#) = 3 ,
 - [nMos36LN](#) = 6 ,
 - [nMos64LN](#) = 7 }
- enum class [DeviceEnumNet](#) {
 - [MCS_DEVICE_ANY](#) ,
 - [MCS_GENERIC_DEVELOPMENT_DEVICE](#) ,
 - [MCS_DEVICE_USB](#) ,
 - [MCS_MCCARD_DEVICE](#) ,
 - [MCS_STG_DEVICE](#) ,
 - [MCS_MC_STIMULUS_DEVICE](#) ,
 - [MCS_MEASUSB_DEVICE](#) ,
 - [MCS_MEA_DEVICE](#) ,
 - [MCS_OCTOPOT_DEVICE](#) ,


```

MCS_TERSENS_DEVICE ,
MCS_PGA_DEVICE ,
MCS_PCX_DEVICE ,
MCS_TCX_DEVICE ,
MCS_FCX_DEVICE ,
MCS_RETINA_LED_DEVICE ,
MCS_MEA_SWITCH_DEVICE ,
MCS_MEA_IMPEDANCE_DEVICE ,
MCS_CHANNELTEST_DEVICE ,
MCS_SW2TO64_DEVICE ,
MCS_RETINA_AMS_DONGLE ,
MCS_PATHIDENT_DEVICE ,
MCS_ROBO_DEVICE ,
MCS_ROBOOCYTE2_DEVICE ,
MCS_ROBOINJECT_DEVICE ,
MCS_HICLAMP_DEVICE ,
MCS_PATCHSERVER_DEVICE ,
MCS_ENCAPSULATOR_DEVICE ,
MCS_MEASURETABLE_DEVICE ,
MCS_FYI_DEVICE ,
MCS_HLA_DEVICE ,
MCS_PPS_DEVICE ,
MCS_PPS5_DEVICE ,
MCS_OKUVISION_STIMULATOR_DEVICE ,
MCS_NF_GEN_DEVICE ,
MCS_SAFEIS_DEVICE ,
MCS_PERISTALTIC_PUMP_DEVICE ,
MCS_EXTERN_BC_TESTER_DEVICE ,
MCS_EXTERN_D_TESTER_DEVICE ,
MCS_SOFTWARE_DONGLE_DEVICE ,
MCS_MEA_CLEAN_DEVICE ,
MCS_MEA_COAT_DEVICE ,
MCS_SMARTIMPLANT_DEVICE ,
MCS_MBC08_DEVICE ,
MCS_PEDOTER_DEVICE ,
MCS_PPC_DEVICE ,
WARNER_VALVE_CONTROL_DEVICE = 7000 ,
WARNER_USSING_DEVICE ,
HEKA_LIH3_DEVICE = 8000 ,
ALA_VC3_DEVICE = 9990 ,
MCS_DEVICE_USB_CYPRESS = 9991 }

```

Enumerates the group of MCS devices to connect to.

- enum class [VendorIdEnumNet](#) {
[Any](#) = -1 ,
[None](#) = 0 ,
[MCS](#) = MCS_VENDOR_ID ,
[PCI](#) = 0x10E8 ,
[Cypress](#) = CYPRESS_VENDOR_ID ,
[ALA_VC3](#) = ALA_VC3_VENDOR_ID }

Enumerates the group of MCS devices to connect to.

- enum class [ProductIdEnumNet](#) {
[Any](#) = -1 ,
[None](#) = 0 ,
[ALA_VC3](#) = ALA_VC3_VENDOR_ID ,
[Cypress_FX1](#) = CY_FX1_PRODUCT_ID ,
[Cypress_FX2](#) = CY_FX2_PRODUCT_ID ,
[Cypress_FX3](#) = CY_FX3_PRODUCT_ID ,
[MC_Card](#) = MCS_PRODUCT_ID_MC_CARD ,

Campden_Ci4600EphysVideoDataIntegrator = MCS_PRODUCT_ID_CAMPDEN_CI4600EPHYS_VIDEO↵
 _DATA_INTEGRATOR ,
 HekaLIH30 = MCS_PRODUCT_ID_HEKA_LIH30 ,
 HekaEPC10Single = MCS_PRODUCT_ID_HEKA_EPC10_SINGLE ,
 HekaEPC10Double = MCS_PRODUCT_ID_HEKA_EPC10_DOUBLE ,
 HekaEPC10Triple = MCS_PRODUCT_ID_HEKA_EPC10_TRIPLE ,
 HekaEPC10Quadro = MCS_PRODUCT_ID_HEKA_EPC10_QUADRO ,
 HekaLIH406 = MCS_PRODUCT_ID_HEKA_LIH_406 ,
 HekaLIH816 = MCS_PRODUCT_ID_HEKA_LIH_816 ,
 HekaITEV100 = MCS_PRODUCT_ID_HEKA_ITEV_100 ,
 HekaPG610 = MCS_PRODUCT_ID_HEKA_PG_610 ,
 HekaPG611 = MCS_PRODUCT_ID_HEKA_PG_611 ,
 HekaPG612 = MCS_PRODUCT_ID_HEKA_PG_612 ,
 HekaPG618 = MCS_PRODUCT_ID_HEKA_PG_618 ,
 HekaPG690 = MCS_PRODUCT_ID_HEKA_PG_690 ,
 HekaEPCLite = MCS_PRODUCT_ID_HEKA_EPC_Lite ,
 STG = MCS_PRODUCT_ID_STG ,
 Octopot = MCS_PRODUCT_ID_OCTOPOT ,
 Tersens = MCS_PRODUCT_ID_TERSENS ,
 Dotriapot = MCS_PRODUCT_ID_DOTRIAPOT ,
 HLA = MCS_PRODUCT_ID_HLA ,
 STG400x = MCS_PRODUCT_ID_STG400x ,
 STG4002 = MCS_PRODUCT_ID_STG4002 ,
 STG4004 = MCS_PRODUCT_ID_STG4004 ,
 STG4008 = MCS_PRODUCT_ID_STG4008 ,
 STG400x_opto = MCS_PRODUCT_ID_STG400x_OPTO ,
 STG4002_opto = MCS_PRODUCT_ID_STG4002_OPTO ,
 STG4004_opto = MCS_PRODUCT_ID_STG4004_OPTO ,
 STG4008_opto = MCS_PRODUCT_ID_STG4008_OPTO ,
 STG5 = MCS_PRODUCT_ID_STG5 ,
 STG3008_FA = MCS_PRODUCT_ID_STG3008_FA ,
 MultiwellOptoStim = MCS_PRODUCT_ID_MULTIWELLOPTOSTIM ,
 Generic = MCS_PRODUCT_ID_GENERIC ,
 PGA = MCS_PRODUCT_ID_PGA ,
 PCX = MCS_PRODUCT_ID_PCX ,
 TCX = MCS_PRODUCT_ID_TCX ,
 FCX = MCS_PRODUCT_ID_FCX ,
 FCB = MCS_PRODUCT_ID_FCB ,
 TC01 = MCS_PRODUCT_ID_TC01 ,
 TC02 = MCS_PRODUCT_ID_TC02 ,
 Retina_LED = MCS_PRODUCT_ID_RETINA_LED ,
 AMS_Dongle = MCS_PRODUCT_ID_RETINA_AMS_DONGLE ,
 OkuVision_Stimulator = MCS_PRODUCT_ID_OKUVISION_STIMULATOR ,
 ExternBCTester = MCS_PRODUCT_ID_RETINAIMPLANT_EXTERNBCTESTER ,
 Triggerbox_IMS = MCS_PRODUCT_ID_RIAG_TRIGGERBOX_IMS ,
 Triggerbox_AMS = MCS_PRODUCT_ID_RIAG_TRIGGERBOX_AMS ,
 Triggerbox_AMS3 = MCS_PRODUCT_ID_RIAG_TRIGGERBOX_AMS3 ,
 ExternDTester = MCS_PRODUCT_ID_RETINAIMPLANT_EXTERNDTESTER ,
 FunkDongleS = MCS_PRODUCT_ID_RIAG_FUNKDONGLES ,
 ExternSTester = MCS_PRODUCT_ID_RIAG_EXTERNSTESTER ,
 DongleS = MCS_PRODUCT_ID_RIAG_DONGLES ,
 Triggerbox_R5 = MCS_PRODUCT_ID_RIAG_TRIGGERBOX_R5 ,
 MEA_Switch = MCS_PRODUCT_ID_MEA_SWITCH ,
 MEA_Impedance = MCS_PRODUCT_ID_MEA_IMPEDANCE ,
 ChannelTest = MCS_PRODUCT_ID_CHANNELTEST ,
 Sw2to64 = MCS_PRODUCT_ID_SW2TO64 ,
 PeristalticPump = MCS_PRODUCT_ID_PERISTALTIC_PUMP ,
 MEA_Switch_2_1 = MCS_PRODUCT_ID_MEA_SWITCH_2_1 ,

```
MEA_Switch_4_2 = MCS_PRODUCT_ID_MEA_SWITCH_4_2 ,
PPS4plus1 = MCS_PRODUCT_ID_PPS4plus1 ,
PPS5 = MCS_PRODUCT_ID_PPS5 ,
PPS2 = MCS_PRODUCT_ID_PPS2 ,
PPS5_DIG = MCS_PRODUCT_ID_PPS5_DIG ,
MEA_Clean = MCS_PRODUCT_ID_MEA_CLEAN ,
MEA_Coat = MCS_PRODUCT_ID_MEA_COAT ,
Multiwell_ICC = MCS_PRODUCT_ID_MULTIWELL_ICC ,
MBC08 = MCS_PRODUCT_ID_MBC08 ,
PPC = MCS_PRODUCT_ID_PPC ,
MEA1060 = MCS_PRODUCT_ID_MEA1060 ,
MEA_Sanofi = MCS_PRODUCT_ID_MEA_SANOFI ,
ME256 = MCS_PRODUCT_ID_ME256 ,
ME128 = MCS_PRODUCT_ID_ME128 ,
ME64 = MCS_PRODUCT_ID_ME64 ,
ME32 = MCS_PRODUCT_ID_ME32 ,
ME16 = MCS_PRODUCT_ID_ME16 ,
MEA2100_Mini_Usb_develop = MCS_PRODUCT_ID_MEA2100_MINI_USB_DEVELOP ,
MEA256 = MCS_PRODUCT_ID_MEA256 ,
MEA2100 = MCS_PRODUCT_ID_MEA2100 ,
MEA2100_32 = MCS_PRODUCT_ID_MEA2100_32 ,
MEA2100_Lite = MCS_PRODUCT_ID_MEA21_LITE ,
Multiwell = MCS_PRODUCT_ID_MULTIWELL ,
MEA2100_256 = MCS_PRODUCT_ID_MEA2100_256 ,
ME2100 = MCS_PRODUCT_ID_ME2100 ,
MEA2100BetaScreen = MCS_PRODUCT_ID_MEA2100_BETA_SCREEN ,
MEA2100_Mini = MCS_PRODUCT_ID_MEA2100_MINI ,
TBSI_Dacq = MCS_PRODUCT_ID_TBSI_DACQ ,
Multiwell_MEA_Mini = MCS_PRODUCT_ID_MULTIWELL_MEA_MINI ,
Whole_Cell_Patch = MCS_PRODUCT_ID_WHOLE_CELL_PATCH ,
eCube = MCS_PRODUCT_ID_ECUBE ,
Graphene_ASIC = MCS_PRODUCT_ID_GRAPHENE_ASIC ,
GE2100 = MCS_PRODUCT_ID_GE2100 ,
Multiboot = MCS_PRODUCT_ID_MULTIBOOT ,
WPA8 = MCS_PRODUCT_ID_WPA8 ,
WPA4 = MCS_PRODUCT_ID_WPA4 ,
WPA16 = MCS_PRODUCT_ID_WPA16 ,
WPA32 = MCS_PRODUCT_ID_WPA32 ,
W2100 = MCS_PRODUCT_ID_W2100 ,
NeuroChip = MCS_PRODUCT_ID_NEUROCHIP ,
UsbTest = MCS_PRODUCT_ID_USB_TEST ,
SoftwareDongle = MCS_PRODUCT_ID_SOFTWAREDONGLE ,
PathIdent = MCS_PRODUCT_ID_PATHIDENT ,
NF_Gen = MCS_PRODUCT_ID_NF_GEN ,
SafeIS = MCS_PRODUCT_ID_SAFEIS ,
Encapsulator = MCS_PRODUCT_ID_ENCAPSULATOR ,
NeurochipConfig = MCS_PRODUCT_ID_NEUROCHIP_CONFIG ,
MeasureTable = MCS_PRODUCT_ID_MEASURETABLE ,
Roboocyte2 = MCS_PRODUCT_ID_ROBOOCYTE2 ,
RoboInject = MCS_PRODUCT_ID_ROBOINJECT ,
HiClamp = MCS_PRODUCT_ID_HICLAMP ,
PatchServer = MCS_PRODUCT_ID_PATCHSERVER ,
Dilutor = MCS_PRODUCT_ID_DILUTOR ,
HiClamp4Uart = MCS_PRODUCT_ID_HICLAMP4UART ,
IM16S16KRA = MCS_PRODUCT_ID_IM16S16KRA ,
IM64KRB = MCS_PRODUCT_ID_IM64KRB ,
IS32KRA = MCS_PRODUCT_ID_IS32KRA ,
IM64KRC = MCS_PRODUCT_ID_IM64KRC ,
```

```

IM16S8KRA = MCS_PRODUCT_ID_IM16S8KRA ,
IM16KRC = MCS_PRODUCT_ID_IM16KRC ,
SmartImplant = MCS_PRODUCT_ID_SMARTIMPLANT ,
PositionImp = MCS_PRODUCT_ID_POSITION_IMP ,
PositionBase = MCS_PRODUCT_ID_POSITION_BASE ,
PositionIICentralUnit = MCS_PRODUCT_ID_POSITIONII_CENTRAL_UNIT ,
PositionIIBase = MCS_PRODUCT_ID_POSITIONII_BASE ,
GrapheneProjectTestDevice = MCS_PRODUCT_ID_GRAPHENE_PROJECT_TEST_DEVICE ,
Pos900 = MCS_PRODUCT_ID_POS900 ,
Neptun = MCS_PRODUCT_ID_NEPTUN ,
Warner_Valve_Control = MCS_PRODUCT_ID_WARNER_VALVE_CONTROL ,
Warner_TEER_Machine = MCS_PRODUCT_ID_WARNER_TEER_MACHINE ,
Warner_Ussing = MCS_PRODUCT_ID_WARNER_USSING }

```

Enumerates the group of MCS devices to connect to.

- enum class `McsBusTypeEnumNet` {
`MCS_ANY_BUS` = -1 ,
`MCS_UNDEFINED_BUS` = 0 ,
`MCS_USB_BUS` ,
`MCS_PCI_BUS` }

Enumerates the bus to use, either USB, PCI or any

- enum class `McsUsbSpeedEnumNet` {
`LowSpeed` = 0 ,
`FullSpeed` = 1 ,
`HighSpeed` = 2 ,
`SuperSpeed` = 3 ,
`UnknownSpeed` = 0xff }

Enumerates the current connection speed of the device

- enum class `CFirmwareDestinationNet` {
`FPGA_NORMAL` = 0 ,
`DSP` = MCSUSB_DEST_DSP ,
`USB` = MCSUSB_DEST_USB ,
`MCU1` = MCSUSB_DEST_MCU1 ,
`MCSBUS1` = MCSUSB_DEST_MCSBUS1 ,
`MCSBUS2` = MCSUSB_DEST_MCSBUS2 ,
`MCSBUS3` = MCSUSB_DEST_MCSBUS3 ,
`MCSBUS4` = MCSUSB_DEST_MCSBUS4 ,
`MCSBUS5` = MCSUSB_DEST_MCSBUS5 ,
`MCSBUS6` = MCSUSB_DEST_MCSBUS6 ,
`MCSBUS7` = MCSUSB_DEST_MCSBUS7 ,
`MCSBUS8` = MCSUSB_DEST_MCSBUS8 ,
`MCSBUS9` = MCSUSB_DEST_MCSBUS9 ,
`MCSBUS10` = MCSUSB_DEST_MCSBUS10 ,
`MCSBUS11` = MCSUSB_DEST_MCSBUS11 ,
`MCSBUS12` = MCSUSB_DEST_MCSBUS12 ,
`MCSBUS13` = MCSUSB_DEST_MCSBUS13 ,
`MCSBUS14` = MCSUSB_DEST_MCSBUS14 ,
`MCSBUS15` = MCSUSB_DEST_MCSBUS15 ,
`MCSBUS0` = MCSUSB_DEST_MCSBUS0 ,
`BUSNUMBER1` = MCSUSB_DEST_BUSNUMBER1 ,
`BUS1MCSBUS1` = MCSUSB_DEST_BUS1_MCSBUS1 ,
`BUS1MCSBUS2` = MCSUSB_DEST_BUS1_MCSBUS2 ,
`BUSNUMBER2` = MCSUSB_DEST_BUSNUMBER2 ,
`BUS2MCSBUS1` = MCSUSB_DEST_BUS2_MCSBUS1 ,
`BUS2MCSBUS2` = MCSUSB_DEST_BUS2_MCSBUS2 ,
`PIC` = MCSUSB_DEST_PIC ,
`PIC2` = MCSUSB_DEST_PIC2 ,
`PIC3` = MCSUSB_DEST_PIC3 ,
`PIC4` = MCSUSB_DEST_PIC4 ,

```

PIC5 = MCSUSB_DEST_PIC5 ,
PIC6 = MCSUSB_DEST_PIC6 ,
PIC7 = MCSUSB_DEST_PIC7 ,
PIC8 = MCSUSB_DEST_PIC8 ,
PIC9 = MCSUSB_DEST_PIC9 ,
PIC10 = MCSUSB_DEST_PIC10 ,
PIC11 = MCSUSB_DEST_PIC11 ,
PIC12 = MCSUSB_DEST_PIC12 ,
ChannelPIC = MCSUSB_DEST_CHANNELPIC ,
Bootstrap = MCSUSB_DEST_BOOTSTRAP ,
BootstrapOtherCypress = MCSUSB_DEST_BOOTSTAP_OTHER_CYPRESS ,
ALTERA = MCSUSB_DEST_ALTERA ,
FPGA2 = MCSUSB_DEST_FPGA2 ,
FPGA3 = MCSUSB_DEST_FPGA3 ,
FPGA4 = MCSUSB_DEST_FPGA4 ,
FPGA5 = MCSUSB_DEST_FPGA5 ,
FPGA6 = MCSUSB_DEST_FPGA6 ,
FPGA7 = MCSUSB_DEST_FPGA7 ,
FPGA8 = MCSUSB_DEST_FPGA8 ,
FPGA9 = MCSUSB_DEST_FPGA9 ,
FPGA10 = MCSUSB_DEST_FPGA10 ,
FPGA11 = MCSUSB_DEST_FPGA11 ,
FPGA12 = MCSUSB_DEST_FPGA12 ,
FPGA13 = MCSUSB_DEST_FPGA13 ,
FPGA14 = MCSUSB_DEST_FPGA14 ,
FPGA15 = MCSUSB_DEST_FPGA15 ,
FPGA16 = MCSUSB_DEST_FPGA16 ,
FPGA_GOLD = XILINX_DEST_GOLDEN ,
ALTERA_GOLD = (MCSUSB_DEST_ALTERA | XILINX_DEST_GOLDEN) ,
FPGA2_GOLD = (MCSUSB_DEST_FPGA2 | XILINX_DEST_GOLDEN) ,
FPGA3_GOLD = (MCSUSB_DEST_FPGA3 | XILINX_DEST_GOLDEN) ,
FPGA4_GOLD = (MCSUSB_DEST_FPGA4 | XILINX_DEST_GOLDEN) ,
FPGA5_GOLD = (MCSUSB_DEST_FPGA5 | XILINX_DEST_GOLDEN) ,
FPGA6_GOLD = (MCSUSB_DEST_FPGA6 | XILINX_DEST_GOLDEN) ,
FPGA7_GOLD = (MCSUSB_DEST_FPGA7 | XILINX_DEST_GOLDEN) ,
FPGA8_GOLD = (MCSUSB_DEST_FPGA8 | XILINX_DEST_GOLDEN) ,
FPGA9_GOLD = (MCSUSB_DEST_FPGA9 | XILINX_DEST_GOLDEN) ,
FPGA10_GOLD = (MCSUSB_DEST_FPGA10 | XILINX_DEST_GOLDEN) ,
FPGA11_GOLD = (MCSUSB_DEST_FPGA11 | XILINX_DEST_GOLDEN) ,
FPGA12_GOLD = (MCSUSB_DEST_FPGA12 | XILINX_DEST_GOLDEN) ,
FPGA13_GOLD = (MCSUSB_DEST_FPGA13 | XILINX_DEST_GOLDEN) ,
FPGA14_GOLD = (MCSUSB_DEST_FPGA14 | XILINX_DEST_GOLDEN) ,
FPGA15_GOLD = (MCSUSB_DEST_FPGA15 | XILINX_DEST_GOLDEN) ,
FPGA16_GOLD = (MCSUSB_DEST_FPGA16 | XILINX_DEST_GOLDEN) ,
FPGA_BASE = XILINX_DEST_BASEIMAGE ,
ALTERA_BASE = (MCSUSB_DEST_ALTERA | XILINX_DEST_BASEIMAGE) ,
FPGA2_BASE = (MCSUSB_DEST_FPGA2 | XILINX_DEST_BASEIMAGE) ,
FPGA3_BASE = (MCSUSB_DEST_FPGA3 | XILINX_DEST_BASEIMAGE) ,
FPGA4_BASE = (MCSUSB_DEST_FPGA4 | XILINX_DEST_BASEIMAGE) ,
FPGA5_BASE = (MCSUSB_DEST_FPGA5 | XILINX_DEST_BASEIMAGE) ,
FPGA6_BASE = (MCSUSB_DEST_FPGA6 | XILINX_DEST_BASEIMAGE) ,
FPGA7_BASE = (MCSUSB_DEST_FPGA7 | XILINX_DEST_BASEIMAGE) ,
FPGA8_BASE = (MCSUSB_DEST_FPGA8 | XILINX_DEST_BASEIMAGE) ,
FPGA9_BASE = (MCSUSB_DEST_FPGA9 | XILINX_DEST_BASEIMAGE) ,
FPGA10_BASE = (MCSUSB_DEST_FPGA10 | XILINX_DEST_BASEIMAGE) ,
FPGA11_BASE = (MCSUSB_DEST_FPGA11 | XILINX_DEST_BASEIMAGE) ,
FPGA12_BASE = (MCSUSB_DEST_FPGA12 | XILINX_DEST_BASEIMAGE) ,
FPGA13_BASE = (MCSUSB_DEST_FPGA13 | XILINX_DEST_BASEIMAGE) ,

```

```

FPGA14_BASE = (MCSUSB_DEST_FPGA14 | XILINX_DEST_BASEIMAGE) ,
FPGA15_BASE = (MCSUSB_DEST_FPGA15 | XILINX_DEST_BASEIMAGE) ,
FPGA16_BASE = (MCSUSB_DEST_FPGA16 | XILINX_DEST_BASEIMAGE) ,
FPGA_BOOTSTRAP = XILINX_DEST_BOOTSTRAP ,
ALTERA_BOOTSTRAP = (MCSUSB_DEST_ALTERA | XILINX_DEST_BOOTSTRAP) ,
DEST_TARGET1 = FLASH_DEST_TARGET1 ,
DEST_TARGET2 = FLASH_DEST_TARGET2 ,
DEST_TARGET3 = FLASH_DEST_TARGET3 ,
DEST_TARGET4 = FLASH_DEST_TARGET4 ,
DEST_TARGET5 = FLASH_DEST_TARGET5 ,
DEST_TARGET6 = FLASH_DEST_TARGET6 ,
DEST_TARGET7 = FLASH_DEST_TARGET7 ,
DEST_TARGET8 = FLASH_DEST_TARGET8 ,
DEST_TARGET9 = FLASH_DEST_TARGET9 ,
DEST_TARGET10 = FLASH_DEST_TARGET10 ,
DEST_TARGET11 = FLASH_DEST_TARGET11 ,
DEST_TARGET12 = FLASH_DEST_TARGET12 ,
DEST_TARGET13 = FLASH_DEST_TARGET13 ,
DEST_TARGET14 = FLASH_DEST_TARGET14 ,
DEST_TARGET15 = FLASH_DEST_TARGET15 ,
DEST_TARGET_MASK = FPGA_DEST_TARGET_MASK ,
DEST_FX3_TARGET_MASK = FX3_DEST_TARGET_MASK ,
ALTERA_TARGET1 = (MCSUSB_DEST_ALTERA | FLASH_DEST_TARGET1) ,
ALTERA_TARGET2 = (MCSUSB_DEST_ALTERA | FLASH_DEST_TARGET2) ,
ALTERA_TARGET3 = (MCSUSB_DEST_ALTERA | FLASH_DEST_TARGET3) ,
USB_TARGET1 = (MCSUSB_DEST_USB | FLASH_DEST_TARGET1) ,
USB_TARGET2 = (MCSUSB_DEST_USB | FLASH_DEST_TARGET2) ,
USB_TARGET3 = (MCSUSB_DEST_USB | FLASH_DEST_TARGET3) ,
UnknownDest = MCSUSB_DEST_UNKNOWN }

```

Enumerates the destination processor for the firmware.

- enum class `DigitalTargetEnumNet` {


```

Digout = (MEA_COMMAND << 16) + MEA_MEA21_DIGOUT_SOURCE ,
Digstream = (MEA_COMMAND << 16) + MEA_MEA21_DIGSTREAM_SOURCE ,
DacqTrigger = (MEA_COMMAND << 16) + MEA_MEA21_DACQTRIGGER_SOURCE ,
StgTrigger = (STG200x_COMMAND << 16) + STG200x_TRIGGER_SOURCE ,
StgListModeTrigger = (STG200x_COMMAND << 16) + STG200x_MEA21_LISTMODE_TRIGGERSOURCE
,
DigOutStimulatorStartTrigger = (MEA_COMMAND << 16) + MEA_DIGOUT_STG_START_TRIGGER_↵
SOURCE ,
DigOutStimulatorStopTrigger = (MEA_COMMAND << 16) + MEA_DIGOUT_STG_STOP_TRIGGER_↵
SOURCE ,
DigStreamToReceiver = (MEA_COMMAND << 16) + MEA_DIGSTREAMTORECEIVER_SOURCE }

```

Enumerates the Digital Targets for Digital Sources

- enum class `DigitalSourceEnumNet` {


```

DigitalInOfOutPort = 0 ,
DigitalIn = 16 ,
DigitalPulse = 32 ,
Feedback = 64 ,
AuxIn = 96 ,
Zero = 98 ,
One = 99 ,
HS1Trigger1Status = 100 ,
HS1Trigger2Status = 102 ,
HS1Trigger3Status = 104 ,
HS1Trigger4Status = 106 ,
HS1Trigger5Status = 108 ,
HS1Trigger6Status = 110 ,
HS1Sideband1 = 112 ,

```

```

HS1Sideband2 = 128 ,
HS1Sideband3 = 144 ,
HS1Sideband4 = 160 ,
HS1Sideband5 = 176 ,
HS1Sideband6 = 192 ,
HS2Trigger1Status = 208 ,
HS2Trigger2Status = 210 ,
HS2Trigger3Status = 212 ,
HS2Trigger4Status = 214 ,
HS2Trigger5Status = 216 ,
HS2Trigger6Status = 218 ,
HS2Sideband1 = 220 ,
HS2Sideband2 = 236 ,
HS2Sideband3 = 252 ,
HS2Sideband4 = 268 ,
HS2Sideband5 = 284 ,
HS2Sideband6 = 300 ,
PulseGenerator = 316 ,
DigitalOutStimulator = 320 ,
DigitalData = 336 ,
DeviceRunStatus = 368 ,
LastPosition = 372 }

```

Enumerates the digital source of the MEA2100 device.

- enum class [W2100DigitalSourceEnumNet](#) {


```

DigitalInOfOutPort = 0 ,
DigitalIn = 16 ,
DigitalPulse = 32 ,
Feedback = 64 ,
AuxIn = 96 ,
Zero = 98 ,
One = 99 ,
PulseGenerator = 100 ,
DigDataFromReceiver = 128 ,
DigitalOutStimulator = 192 ,
DigitalData = 208 ,
DeviceRunStatus = 240 ,
DigStreamFromReceiver = 256 ,
LastPosition = 320 }

```

Enumerates the digital source of the W2100 device.

- enum class [SCUDigitalSourceEnumNet](#) {


```

DigitalInOfOutPort = (0x00 << 8) + 0 ,
DigitalIn = (0x00 << 8) + 16 ,
DigitalPulse = (0x01 << 8) ,
Feedback = (0x02 << 8) ,
AuxIn = (0x03 << 8) + 0 ,
Zero = (0x03 << 8) + 2 ,
One = (0x03 << 8) + 3 ,
PulseGenerator = (0x03 << 8) + 8 ,
DigitalOutStimulator = (0x03 << 8) + 16 ,
DigitalData = (0x04 << 8) ,
DeviceRunStatus = (0x05 << 8) + 0 ,
HS1Trigger1Status = (0x40 << 8) + 0 ,
HS1Trigger2Status = (0x40 << 8) + 2 ,
HS1Trigger3Status = (0x40 << 8) + 4 ,
HS1Trigger4Status = (0x40 << 8) + 6 ,
HS1Trigger5Status = (0x40 << 8) + 8 ,
HS1Trigger6Status = (0x40 << 8) + 10 ,
HS1Trigger7Status = (0x40 << 8) + 12 ,

```



```

HS1Trigger8Status = (0x40 << 8) + 14 ,
HS1Trigger9Status = (0x40 << 8) + 16 ,
HS1Trigger10Status = (0x40 << 8) + 18 ,
HS1Trigger11Status = (0x40 << 8) + 20 ,
HS1Trigger12Status = (0x40 << 8) + 22 ,
HS1Sideband1 = (0x42 << 8) ,
HS1Sideband2 = (0x43 << 8) ,
HS1Sideband3 = (0x44 << 8) ,
HS1Sideband4 = (0x45 << 8) ,
HS1Sideband5 = (0x46 << 8) ,
HS1Sideband6 = (0x47 << 8) ,
HS1Sideband7 = (0x48 << 8) ,
HS1Sideband8 = (0x49 << 8) ,
HS1Sideband9 = (0x4A << 8) ,
HS1Sideband10 = (0x4B << 8) ,
HS1Sideband11 = (0x4C << 8) ,
HS1Sideband12 = (0x4D << 8) ,
HS2Trigger1Status = (0x80 << 8) + 0 ,
HS2Trigger2Status = (0x80 << 8) + 2 ,
HS2Trigger3Status = (0x80 << 8) + 4 ,
HS2Trigger4Status = (0x80 << 8) + 6 ,
HS2Trigger5Status = (0x80 << 8) + 8 ,
HS2Trigger6Status = (0x80 << 8) + 10 ,
HS2Trigger7Status = (0x80 << 8) + 12 ,
HS2Trigger8Status = (0x80 << 8) + 14 ,
HS2Trigger9Status = (0x80 << 8) + 16 ,
HS2Trigger10Status = (0x80 << 8) + 18 ,
HS2Trigger11Status = (0x80 << 8) + 20 ,
HS2Trigger12Status = (0x80 << 8) + 22 ,
HS2Sideband1 = (0x82 << 8) ,
HS2Sideband2 = (0x83 << 8) ,
HS2Sideband3 = (0x84 << 8) ,
HS2Sideband4 = (0x85 << 8) ,
HS2Sideband5 = (0x86 << 8) ,
HS2Sideband6 = (0x87 << 8) ,
HS2Sideband7 = (0x88 << 8) ,
HS2Sideband8 = (0x89 << 8) ,
HS2Sideband9 = (0x8A << 8) ,
HS2Sideband10 = (0x8B << 8) ,
HS2Sideband11 = (0x8C << 8) ,
HS2Sideband12 = (0x8D << 8) ,
LastPosition = (0xFF << 8) }

```

Enumerates the digital source of the SCU device.

- enum class MEA2100_256DigitalSourceEnumNet {

```

DigitalInOfOutPort = (0x00 << 8) + 0 ,
DigitalIn = (0x00 << 8) + 16 ,
DigitalPulse = (0x01 << 8) ,
Feedback = (0x02 << 8) ,
AuxIn = (0x03 << 8) + 0 ,
Zero = (0x03 << 8) + 2 ,
One = (0x03 << 8) + 3 ,
DeviceRunStatus = (0x03 << 8) + 4 ,
PulseGenerator = (0x03 << 8) + 8 ,
DigitalOutStimulator = (0x03 << 8) + 16 ,
DigitalData = (0x04 << 8) ,
HS1Trigger1Status = (0x40 << 8) + 0 ,
HS1Trigger2Status = (0x40 << 8) + 2 ,
HS1Trigger3Status = (0x40 << 8) + 4 ,

```



```
HS1Trigger4Status = (0x40 << 8) + 6 ,
HS1Trigger5Status = (0x40 << 8) + 8 ,
HS1Trigger6Status = (0x40 << 8) + 10 ,
HS1Trigger7Status = (0x40 << 8) + 12 ,
HS1Trigger8Status = (0x40 << 8) + 14 ,
HS1Trigger9Status = (0x40 << 8) + 16 ,
HS1Trigger10Status = (0x40 << 8) + 18 ,
HS1Trigger11Status = (0x40 << 8) + 20 ,
HS1Trigger12Status = (0x40 << 8) + 22 ,
HS1Trigger13Status = (0x40 << 8) + 24 ,
HS1Trigger14Status = (0x40 << 8) + 26 ,
HS1Trigger15Status = (0x40 << 8) + 28 ,
HS1Trigger16Status = (0x40 << 8) + 30 ,
HS1Trigger17Status = (0x41 << 8) + 0 ,
HS1Trigger18Status = (0x41 << 8) + 2 ,
HS1Sideband1 = (0x42 << 8) ,
HS1Sideband2 = (0x43 << 8) ,
HS1Sideband3 = (0x44 << 8) ,
HS1Sideband4 = (0x45 << 8) ,
HS1Sideband5 = (0x46 << 8) ,
HS1Sideband6 = (0x47 << 8) ,
HS1Sideband7 = (0x48 << 8) ,
HS1Sideband8 = (0x49 << 8) ,
HS1Sideband9 = (0x4A << 8) ,
HS1Sideband10 = (0x4B << 8) ,
HS1Sideband11 = (0x4C << 8) ,
HS1Sideband12 = (0x4D << 8) ,
HS1Sideband13 = (0x4E << 8) ,
HS1Sideband14 = (0x4F << 8) ,
HS1Sideband15 = (0x50 << 8) ,
HS1Sideband16 = (0x51 << 8) ,
HS1Sideband17 = (0x52 << 8) ,
HS1Sideband18 = (0x53 << 8) ,
HS2Trigger1Status = (0x80 << 8) + 0 ,
HS2Trigger2Status = (0x80 << 8) + 2 ,
HS2Trigger3Status = (0x80 << 8) + 4 ,
HS2Trigger4Status = (0x80 << 8) + 6 ,
HS2Trigger5Status = (0x80 << 8) + 8 ,
HS2Trigger6Status = (0x80 << 8) + 10 ,
HS2Trigger7Status = (0x80 << 8) + 12 ,
HS2Trigger8Status = (0x80 << 8) + 14 ,
HS2Trigger9Status = (0x80 << 8) + 16 ,
HS2Trigger10Status = (0x80 << 8) + 18 ,
HS2Trigger11Status = (0x80 << 8) + 20 ,
HS2Trigger12Status = (0x80 << 8) + 22 ,
HS2Trigger13Status = (0x80 << 8) + 24 ,
HS2Trigger14Status = (0x80 << 8) + 26 ,
HS2Trigger15Status = (0x80 << 8) + 28 ,
HS2Trigger16Status = (0x80 << 8) + 30 ,
HS2Trigger17Status = (0x81 << 8) + 0 ,
HS2Trigger18Status = (0x81 << 8) + 2 ,
HS2Sideband1 = (0x82 << 8) ,
HS2Sideband2 = (0x83 << 8) ,
HS2Sideband3 = (0x84 << 8) ,
HS2Sideband4 = (0x85 << 8) ,
HS2Sideband5 = (0x86 << 8) ,
HS2Sideband6 = (0x87 << 8) ,
HS2Sideband7 = (0x88 << 8) ,
```

```

HS2Sideband8 = (0x89 << 8) ,
HS2Sideband9 = (0x8A << 8) ,
HS2Sideband10 = (0x8B << 8) ,
HS2Sideband11 = (0x8C << 8) ,
HS2Sideband12 = (0x8D << 8) ,
HS2Sideband13 = (0x8E << 8) ,
HS2Sideband14 = (0x8F << 8) ,
HS2Sideband15 = (0x90 << 8) ,
HS2Sideband16 = (0x91 << 8) ,
HS2Sideband17 = (0x92 << 8) ,
HS2Sideband18 = (0x93 << 8) ,
LastPosition = (0xFF << 8) }

```

Enumerates the digital source of the MEA2100-256 device.

- enum class `TBSI_DACQDigitalSourceEnumNet` {
 `DigitalInOfOutPort` = (0x00 << 8) + 0 ,
 `DigitalIn` = (0x00 << 8) + 16 ,
 `DigitalPulse` = (0x01 << 8) ,
 `Feedback` = (0x02 << 8) ,
 `AuxIn` = (0x03 << 8) + 0 ,
 `Zero` = (0x03 << 8) + 2 ,
 `One` = (0x03 << 8) + 3 ,
 `DeviceRunStatus` = (0x03 << 8) + 4 ,
 `PulseGenerator` = (0x03 << 8) + 8 ,
 `DigitalOutStimulator` = (0x03 << 8) + 16 ,
 `DigitalData` = (0x04 << 8) ,
 `HS1DigitalData1` = (0x30 << 8) ,
 `HS2DigitalData1` = (0x70 << 8) ,
 `LastPosition` = (0xFF << 8) }

Enumerates the digital source of the TBSI-DACQ device.

- enum class `TriggerSourceEnumNet` {
 `tsNone` = 0 ,
 `tsDigitalIn1` = 1 ,
 `tsDigitalIn2` = 2 ,
 `tsDigitalIn3` = 3 ,
 `tsDigitalIn4` = 4 ,
 `tsDigitalIn5` = 5 ,
 `tsDigitalIn6` = 6 ,
 `tsDigitalIn7` = 7 ,
 `tsDigitalIn8` = 8 ,
 `tsDigitalIn9` = 9 ,
 `tsDigitalIn10` = 10 ,
 `tsDigitalIn11` = 11 ,
 `tsDigitalIn12` = 12 ,
 `tsDigitalIn13` = 13 ,
 `tsDigitalIn14` = 14 ,
 `tsDigitalIn15` = 15 ,
 `tsDigitalIn16` = 16 ,
 `tsDigitalIn17` = 17 ,
 `tsDigitalIn18` = 18 ,
 `tsDigitalIn19` = 19 ,
 `tsDigitalIn20` = 20 ,
 `tsDigitalIn21` = 21 ,
 `tsDigitalIn22` = 22 ,
 `tsDigitalIn23` = 23 ,
 `tsDigitalIn24` = 24 ,
 `tsDigitalIn25` = 25 ,
 `tsDigitalIn26` = 26 ,
 `tsDigitalIn27` = 27 ,
 }

```
tsDigitalIn28 = 28 ,  
tsDigitalIn29 = 29 ,  
tsDigitalIn30 = 30 ,  
tsDigitalIn31 = 31 ,  
tsDigitalIn32 = 32 ,  
tsFeedback1 = 33 ,  
tsFeedback2 = 34 ,  
tsFeedback3 = 35 ,  
tsFeedback4 = 36 ,  
tsFeedback5 = 37 ,  
tsFeedback6 = 38 ,  
tsFeedback7 = 39 ,  
tsFeedback8 = 40 ,  
tsFeedback9 = 41 ,  
tsFeedback10 = 42 ,  
tsFeedback11 = 43 ,  
tsFeedback12 = 44 ,  
tsFeedback13 = 45 ,  
tsFeedback14 = 46 ,  
tsFeedback15 = 47 ,  
tsFeedback16 = 48 ,  
tsFeedback17 = 49 ,  
tsFeedback18 = 50 ,  
tsFeedback19 = 51 ,  
tsFeedback20 = 52 ,  
tsFeedback21 = 53 ,  
tsFeedback22 = 54 ,  
tsFeedback23 = 55 ,  
tsFeedback24 = 56 ,  
tsFeedback25 = 57 ,  
tsFeedback26 = 58 ,  
tsFeedback27 = 59 ,  
tsFeedback28 = 60 ,  
tsFeedback29 = 61 ,  
tsFeedback30 = 62 ,  
tsFeedback31 = 63 ,  
tsFeedback32 = 64 ,  
tsAuxIn1 = 65 ,  
tsAuxIn2 = 66 ,  
tsDigitalPuse0 = 67 ,  
tsDigitalPuse1 = 68 ,  
tsDigitalPuse2 = 69 ,  
tsDigitalPuse3 = 70 ,  
tsDigitalPuse4 = 71 ,  
tsDigitalPuse5 = 72 ,  
tsDigitalPuse6 = 73 ,  
tsDigitalPuse7 = 74 ,  
tsDigitalPuse8 = 75 ,  
tsDigitalPuse9 = 76 ,  
tsDigitalPuse10 = 77 ,  
tsDigitalPuse11 = 78 ,  
tsDigitalPuse12 = 79 ,  
tsDigitalPuse13 = 80 ,  
tsDigitalPuse14 = 81 ,  
tsDigitalPuse15 = 82 ,  
tsDigitalPuse16 = 83 ,  
tsDigitalPuse17 = 84 ,  
tsDigitalPuse18 = 85 ,
```

```

tsDigitalPuse19 = 86 ,
tsDigitalPuse20 = 87 ,
tsDigitalPuse21 = 88 ,
tsDigitalPuse22 = 89 ,
tsDigitalPuse23 = 90 ,
tsDigitalPuse24 = 91 ,
tsDigitalPuse25 = 92 ,
tsDigitalPuse26 = 93 ,
tsDigitalPuse27 = 94 ,
tsDigitalPuse28 = 95 ,
tsDigitalPuse29 = 96 ,
tsDigitalPuse30 = 97 ,
tsDigitalPuse31 = 98 ,
tsTriggered = 99 ,
tsSidebandBit8 = 100 ,
tsDACQCy1Dev1Runs = 101 ,
tsDACQCy1Dev2Runs = 102 ,
tsDACQCy2Dev1Runs = 103 ,
tsDACQCy2Dev2Runs = 104 }

```

Enumerates the trigger source of the MEA2100 device.

- enum class [AnalogSourceEnumNet](#) {
[AnalogSource_HS1](#) ,
[AnalogSource_HS2](#) ,
[AnalogSource_IF](#) }

Enumerates the analog source of the MEA2100 device.

- enum class [Stg200xTriggerStatusEnumNet](#) {
[Idle](#) = 0 ,
[Running](#) = 1 ,
[Finished](#) = 2 ,
[Armed](#) = 3 }

Enumerates the STG download mode trigger status

- enum class [RetriggerActionEnumNet](#) {
[raStop](#) = STG200x_RETRIGGER_STOP ,
[raRestart](#) = STG200x_RETRIGGER_RESTART ,
[raIgnore](#) = STG200x_RETRIGGER_IGNORE ,
[raGate](#) = STG200x_RETRIGGER_GATEMODE ,
[raSingle](#) = STG200x_RETRIGGER_SINGLE }

Enumerates possible retrigger actions for STG200x devices.

- enum class [Stg200xSegmentFlagsEnumNet](#) {
[None](#) = 0 ,
[UpdateTrigger](#) = SEGMENTFLAGS_UPDATETRIGGER ,
[DownloadOnly](#) = SEGMENTFLAGS_DOWNLOADONLY ,
[TriggerOnly](#) = SEGMENTFLAGS_TRIGGERONLY ,
[SyncStart](#) = SEGMENTFLAGS_SYNCSTART }

Enumerates Segmentflag options for STG400x devices.

- enum class [Stg200xDigoutModeEnumNet](#) {
[Monitor](#) = STG200x_DIGOUTMODE_MONITOR ,
[Manual](#) = STG200x_DIGOUTMODE_MANUAL ,
[SYNCOUT1](#) = STG200x_DIGOUTMODE_SYNCOUT1 ,
[SYNCOUT2](#) = STG200x_DIGOUTMODE_SYNCOUT2 ,
[SYNCOUT3](#) = STG200x_DIGOUTMODE_SYNCOUT3 ,
[SYNCOUT4](#) = STG200x_DIGOUTMODE_SYNCOUT4 ,
[SYNCOUT5](#) = STG200x_DIGOUTMODE_SYNCOUT5 ,
[SYNCOUT6](#) = STG200x_DIGOUTMODE_SYNCOUT6 ,
[SYNCOUT7](#) = STG200x_DIGOUTMODE_SYNCOUT7 ,
[SYNCOUT8](#) = STG200x_DIGOUTMODE_SYNCOUT8 }

Enumerates the DigoutMode on STG400x devices.

- enum class [DigitalStimulatorTriggerSlopeEnumNet](#) {
[Falling](#) = 0 ,
[Rising](#) = 1 }
Enumerates start/stop conditions for DigOut/DigStim trigger. /summary>
- enum class [DigitalStimulatorTriggerEventEnumNet](#) {
[Start](#) = 0 ,
[Stop](#) = 1 }
Enumerates start/stop event for DigOut/DigStim trigger. /summary>
- enum class [AdapterTypeEnumNet](#) {
[None](#) = 0 ,
[MEA60](#) = 1 ,
[MEA2x60](#) = 2 ,
[MEA120](#) = 3 ,
[MEA32](#) = 4 ,
[MEA2x32](#) = 5 ,
[Multiwell96](#) = 6 ,
[WirelessTestAdapter](#) = 7 ,
[MEA252](#) = 8 ,
[MEA_2_252_2](#) = 9 ,
[MEA_2_252_2_6Well](#) = 10 ,
[MEA_2_252_2_9Well](#) = 11 ,
[MEA_2_252_2_Test](#) = 12 ,
[TBSI_5](#) = 13 ,
[TBSI_15](#) = 14 ,
[TBSI_31](#) = 15 ,
[TBSI_63](#) = 16 ,
[TBSI_127](#) = 17 ,
[TBSI_Reserved](#) = 18 ,
[Ci4600Intan](#) = 20 ,
[Unknown](#) = ADAPTER_TYPE_UNKOWN ,
[NotApplicable](#) = ADAPTER_TYPE_ENUM_NOT_APPLICABLE }
Enumerates the adapter type of the MEA2100 device.
- enum class [MeaLayoutEnumNet](#) {
[mlUnknown](#) = 0 ,
[mlMEA60](#) = 1 }
Enumerates the MEA layout of the MEA2100 device.
- enum class [DataModeEnumNet](#) {
[Unsigned_16bit](#) = 0 ,
[Unsigned_24bit](#) = 2 ,
[Unsigned_32bit](#) = 3 ,
[Signed_16bit](#) = 8 ,
[Signed_24bit](#) = 10 ,
[Signed_32bit](#) = 11 }
Enumerates the data mode of the device, either 16, 24 or 32 bit, can be signed or unsigned.
- enum class [SampleSizeNet](#) {
[SampleSize16Unsigned](#) = 2 ,
[SampleSize16Signed](#) = 2 + 0x100 ,
[SampleSize24Unsigned](#) = 3 ,
[SampleSize24Signed](#) = 3 + 0x100 ,
[SampleSize32Unsigned](#) = 4 ,
[SampleSize32Signed](#) = 4 + 0x100 ,
[SampleSize64Unsigned](#) = 8 ,
[SampleSize64Signed](#) = 8 + 0x100 }
Enumerates the data format for ChannelBlock functions.
- enum class [SampleDstSizeNet](#) {
[SampleDstSize16](#) = 2 ,
[SampleDstSize32](#) = 4 }

Enumerates the destination data format for ChannelBlock functions.

- enum class `TcxDeviceTypeEnumNet` {
`Unknown` = 0 ,
`Regular` = 1 ,
`BMI` = 2 ,
`Nanion` = 3 ,
`Warner` = 4 }

Enumerates the type of TCX devices.

- enum class `TcxSensorTypeEnumNet` {
`Reserved5` = 0 ,
`Reserved4` = 1 ,
`Reserved3` = 2 ,
`Reserved2` = 3 ,
`Reserved1` = 4 ,
`NTC10K` = 5 ,
`PT1000` = 6 ,
`PT100` = 7 }

Enumerates the sensor types for TCX devices

- enum class `STG_DestinationEnumNet` {
`channeldata_voltage` ,
`channeldata_current` ,
`syncoutdata` ,
`channeldata_positive_voltage` ,
`channeldata_positive_current` ,
`rawdata` ,
`channeldata_current_own_sync` ,
`channeldata_positive_current_own_sync` ,
`channeldata_current_own_boost_gnd_sync` ,
`channeldata_positive_current_own_boost_gnd_sync` ,
`channeldata_current_always_boost` ,
`channeldata_current_always_boost_own_sync` }

Enumerates the destination for STG downloads.

- enum class `ElectrodeModeEnumNet` {
`emAutomatic` = 0 ,
`emManual` = 3 }

Enumerates the mode of each electrode, can be automatic or manual. In automatic mode, the blanking of the electrode is controlled by the sideband signal, in manual mode, the stimulation configuration is independant of the sideband signal.

- enum class `ElectrodeDacMuxEnumNet` {
`Ground` = 0 ,
`Stg1` = 1 ,
`Stg2` = 2 ,
`Stg3` = 3 }

Enumerates the setting of the Stimulation DAC Multiplexer.

- enum class `DacqGroupChannelEnumNet` {
`HeadstageElectrodeGroup` = 0x00 ,
`InterfaceADCGroup` = INTERFACEANALOGCHANNELSGROUP ,
`DSPDataGroup` = DSPDATACHANNELSGROUP ,
`Headstage1NCBathCurrentGroup` = 0x30 ,
`Headstage1NCCol2CurrentGroup` = 0x31 ,
`Headstage1NChipTempGroup` = 0x32 ,
`STG1DACSignalGroup` = 0x38 ,
`LIH30UserADCGroup` = 0x50 ,
`LIH30TestADCGroup` = 0x51 ,
`LIH30ADCModulesGroup` = 0x52 ,
`IFDigChannelsGroup` = INTERFACEDIGITALCHANNELSGROUP ,
`STG1SidebandsGroup` = 0x90 ,

```

STG1TriggerStatusGroup = 0x91 ,
DACQ1DigitalGroup = 0xA0 ,
AudioTestChannelGroup = AUDIOTESTCHANNELGROUP ,
PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }

```

Enumerates the Channel Groups of Datastream

- enum class [W2100DacqGroupChannelEnumNet](#) {
 - [InterfaceADCGroup](#) = INTERFACEANALOGCHANNELSGROUP ,
 - [DSPDataGroup](#) = DSPDATACHANNELSGROUP ,
 - [WirelessHeadStageAnalogRE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 ,
 - [WirelessHeadStageStatusRE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 ,
 - [WirelessHeadStageAnalogRE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 ,
 - [WirelessHeadStageStatusRE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 ,
 - [WirelessHeadStageAnalogRE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 ,
 - [WirelessHeadStageStatusRE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 ,
 - [WirelessHeadStageAnalogRE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 ,
 - [WirelessHeadStageStatusRE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 ,
 - [WirelessHeadStageAnalogRE2HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 ,
 - [WirelessHeadStageStatusRE2HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 ,
 - [WirelessHeadStageAnalogRE2HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0 ,
 - [WirelessHeadStageStatusRE2HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 ,
 - [WirelessHeadStageAnalogRE2HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0 ,
 - [WirelessHeadStageStatusRE2HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 ,
 - [WirelessHeadStageAnalogRE2HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0 ,
 - [WirelessHeadStageStatusRE2HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 ,
 - [WirelessHeadStageGyroDataRE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE2HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE2HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE2HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE2HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE2HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE2HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 16 ,
 - [WirelessHeadStageGyroDataRE2HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0 + 16 ,
 - [WirelessHeadStageAccDataRE2HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 16 ,
 - [WirelessHeadStageOptoStimCurrentRE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 + 32 ,
 - [WirelessHeadStageReservedARE1HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 32 ,
 - [WirelessHeadStageOptoStimCurrentRE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 + 32 ,
 - [WirelessHeadStageReservedARE1HS2](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 32 ,
 - [WirelessHeadStageOptoStimCurrentRE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 + 32 ,
 - [WirelessHeadStageReservedARE1HS3](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 32 ,
 - [WirelessHeadStageOptoStimCurrentRE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 + 32 ,
 - [WirelessHeadStageReservedARE1HS4](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 32 ,
 - [WirelessHeadStageOptoStimCurrentRE2HS1](#) = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 + 32 ,

```

WirelessHeadStageReservedARE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 32
,
WirelessHeadStageOptoStimCurrentRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 +
0 + 32 ,
WirelessHeadStageReservedARE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 32
,
WirelessHeadStageOptoStimCurrentRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 +
0 + 32 ,
WirelessHeadStageReservedARE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 32
,
WirelessHeadStageOptoStimCurrentRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 +
0 + 32 ,
WirelessHeadStageReservedARE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 32
,
WirelessHeadStageReservedBRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 + 48
,
WirelessHeadStageReservedCRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 48
,
WirelessHeadStageReservedBRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 + 48
,
WirelessHeadStageReservedCRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 48
,
WirelessHeadStageReservedBRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 + 48
,
WirelessHeadStageReservedCRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 48
,
WirelessHeadStageReservedBRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 + 48
,
WirelessHeadStageReservedCRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 48
,
WirelessHeadStageReservedBRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 + 48
,
WirelessHeadStageReservedCRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 48
,
WirelessHeadStageReservedBRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0 + 48
,
WirelessHeadStageReservedCRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 48
,
WirelessHeadStageReservedBRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0 + 48
,
WirelessHeadStageReservedCRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 48
,
WirelessHeadStageReservedBRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0 + 48
,
WirelessHeadStageReservedCRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 48
,
IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP ,
AudioTestChannelGroup = AUDIOTESTCHANNELGROUP ,
PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }

```

Enumerates the W2100 Channel Groups of Datastream

- enum class SCUDacqGroupChannelEnumNet {
 SCU1ElectrodeGroupHS1 = 0x00 ,
 SCU1ElectrodeGroupHS2 = 0x01 ,
 SCU1ElectrodeGroupHS3 = 0x02 ,
 SCU1ElectrodeGroupHS4 = 0x03 ,
 SCU2ElectrodeGroupHS1 = 0x08 ,
 SCU2ElectrodeGroupHS2 = 0x09 ,
 SCU2ElectrodeGroupHS3 = 0x0A ,


```

SCU2ElectrodeGroupHS4 = 0x0B ,
InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP ,
STG1DACSignalGroup = 0x40 ,
STG2DACSignalGroup = 0x41 ,
DSPAnalogGroup = DSPDATACHANNELSGROUP ,
DSPDigitalGroup = 0xA0 ,
IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP ,
STG1TriggerStatusGroup = 0x90 ,
STG1SidebandsGroup = 0x91 ,
STG2TriggerStatusGroup = 0x98 ,
STG2SidebandsGroup = 0x99 ,
AudioTestChannelGroup = AUDIOTESTCHANNELGROUP ,
PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }

```

Enumerates the SCU Channel Groups of Datastream

- enum class **MEA2100_256DacqGroupChannelEnumNet** {
 HS1ElectrodeGroup = 0x00 ,
 HS2ElectrodeGroup = 0x08 ,
 InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP ,
 STG1DACSignalGroup = 0x40 ,
 STG2DACSignalGroup = 0x41 ,
 DSPAnalogGroup = DSPDATACHANNELSGROUP ,
 DSPDigitalGroup = 0xA0 ,
 IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP ,
 STG1TriggerStatusGroup = 0x90 ,
 STG1SidebandsGroup = 0x91 ,
 STG2TriggerStatusGroup = 0x98 ,
 STG2SidebandsGroup = 0x99 ,
 AudioTestChannelGroup = AUDIOTESTCHANNELGROUP ,
 PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }

Enumerates the MEA2100-256 Channel Groups of Datastream

- enum class **DacqMeaGroupTypeEnumNet** {
 AnalogGroup = ANALOG_GROUP ,
 DigitalGroup = DIGITAL_GROUP ,
 FrameContextGroup = FRAME_CONTEXT_GROUP }

Enumerations of CMOS MEA Groups to detect wether it is an Analog, Digital or Frame Context Group

- enum class **CMOSMeaValueUnitEnumNet** {
 NoUnit = 0x00 ,
 NanoVolt = 0x11 ,
 PicoAmpere = 0x21 ,
 NanoAmpere = 0x22 ,
 MicroAmpere = 0x23 ,
 MilliDegreeCelsius = 0x31 }

Enumerations of CMOS MEA Units of Values in Data stream

- enum class **CMOSMeaInterfaceADCEnumNet** {
 IFChannel1 = 0x01 ,
 IFChannel2 = 0x02 ,
 IFChannel3 = 0x04 ,
 IFChannel4 = 0x08 ,
 IFChannel5 = 0x10 ,
 IFChannel6 = 0x20 ,
 IFChannel7 = 0x40 ,
 IFChannel8 = 0x80 }

Enumerations of CMOS MEA IF Analog Channels Group Bitmask

- enum class **CMOSMeaHeadstage1NCBathCurrentEnumNet** { **NCBathCurrent** = 0x01 }

Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

- enum class **CMOSMeaHeadstage1NCCol2CurrentEnumNet** { **NCCol2Current** = 0x01 }

Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

- enum class `CMOSMeaHeadstage1NChipTempEnumNet` { `NChipTemperature` = 0x01 }

Enumerations of CMOS MEA HS Temperature Monitoring Channels Group Bitmask

- enum class `CMOSMeaSTG1DACSignalEnumNet` {
`DAC1Channel` = 0x01 ,
`DAC2Channel` = 0x02 ,
`DAC3Channel` = 0x04 ,
`DAC4Channel` = 0x08 }

Enumerations of CMOS MEA DAC Stimulation Channels Group Bitmask

- enum class `CMOSMeaIFDigChannelEnumNet` {
`DigitalMux` = 0x01 ,
`DigitalInPort` = 0x02 ,
`DigitalOutReg` = 0x04 ,
`FeedbackReg` = 0x08 ,
`DigitalReg` = 0x10 ,
`AuxPort` = 0x20 }

Enumerations of CMOS MEA IF Digital Channels Group Bitmask

- enum class `CMOSMeaHS1SidebandEnumNet` {
`SBSVector1` = 0x01 ,
`SBSVector2` = 0x02 ,
`SBSVector3` = 0x04 ,
`SBSVector4` = 0x08 }

Enumerations of CMOS MEA HS STG Sideband Channels Group Bitmask

- enum class `CMOSMeaHS1TriggerStatusEnumNet` {
`TriggerStatus1` = 0x01 ,
`TriggerStatus2` = 0x02 ,
`TriggerStatus3` = 0x04 ,
`TriggerStatus4` = 0x08 }

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

- enum class `AnalogUnitEnumNet` {
`Unknown` ,
`Volt` ,
`Ampere` ,
`Kelvin` }
- enum class `CMOSMeaPacketFrameContextGroupEnumNet` {
`SOFAndCTRLword` = 0x01 ,
`ChecksumAndPacketCounter` = 0x02 ,
`Timestamp` = 0x04 ,
`EOFAndCRC` = 0x08 }

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

- enum class `CMOSMeaBathModeEnumNet` {
`Ground` = 0x02 ,
`Stimulation` = 0x01 ,
`CurrentMeasure` = 0x00 }

Enumerations of CMOS MEA Bath Mode

- enum class `PatchServAdcModeEnumNet` {
`Normal` = 0 ,
`CatchAmp` = 1 }
- enum class `RoboCurrentModeEnumNet` {
`Off` = ROBO_CURRENT_OFFMODE ,
`Break` = ROBO_CURRENT_BREAKMODE ,
`Standby` = ROBO_CURRENT_STANDBYMODE ,
`Reference` = ROBO_CURRENT_REFERENCEMODE ,
`Movement` = ROBO_CURRENT_MOVEMENTMODE }

- enum class [TeerClampModeEnumNet](#) {
[ClampModeVoltage](#) = 0 ,
[ClampModeCurrent](#) = 1 ,
[ClampModeOpen](#) = 2 ,
[ClampModeInternalCalibration](#) = 3 }
- enum class [TeerWaveformEnumNet](#) {
[Rectangle](#) = 0 ,
[Sine](#) = 1 }
- enum class [UssingClampModeEnumNet](#) {
[VoltageClamp](#) = 1 ,
[CurrentClamp](#) = 2 ,
[OpenClamp](#) = 3 ,
[Standby](#) = 4 ,
[ElectrodeOffset](#) = 5 }
- enum class [UssingUnitEnumNet](#) {
[Volt](#) = 0 ,
[Ampere](#) = 1 ,
[State](#) = 2 }
- enum class [PlateClampEnumNet](#) {
[Close](#) = 0 ,
[Open](#) = 1 ,
[Stop](#) = 2 }
- enum class [PlateClampLockEnumNet](#) {
[Lock](#) = 0 ,
[Unlock](#) = 1 }
- enum class [MultiwellPlateTypeEnumNet](#) {
[Plate_Dummy](#) = HS_PLATETYPE_0 ,
[Plate_24W700_100FMA](#) = 1 ,
[Plate_24W030MGA](#) = 2 ,
[Plate_72W500_100PMA](#) = 3 ,
[Plate_72W500_100FMA](#) = 5 ,
[Plate_24W700_100FMB](#) = HS_PLATETYPE_6 ,
[Plate_96W700_100FMA](#) = HS_PLATETYPE_7 ,
[Plate_96W300_80_1152FMA](#) = HS_PLATETYPE_33 ,
[Plate_96W400_80_1152FMB](#) = HS_PLATETYPE_36 ,
[Plate_24W300_30_1152GBA](#) = HS_PLATETYPE_40 ,
[Plate_24W700_100FMC](#) = HS_PLATETYPE_44 ,
[Plate_96W700_100FMB](#) = HS_PLATETYPE_48 ,
[Plate_96W700_100GBC](#) = HS_PLATETYPE_49 ,
[Plate_96W700_100GBD](#) = HS_PLATETYPE_51 ,
[Plate_24W700_100PBA](#) = HS_PLATETYPE_60 ,
[Plate_Dummy_126](#) = HS_PLATETYPE_126 ,
[Plate_24W300_30GMA](#) = HS_PLATETYPE_193 ,
[Plate_96W700_100GMA](#) = 194 ,
[Plate_24W300_30GBA](#) = HS_PLATETYPE_195 ,
[Plate_96W700_100GBA](#) = HS_PLATETYPE_224 ,
[Plate_24W300_30GBB](#) = HS_PLATETYPE_232 ,
[Plate_96W700_100GBB](#) = HS_PLATETYPE_244 ,
[No_Plate](#) = 255 }
- enum class [FpgaldEnumNet](#) {
[DeviceNotConnected](#) = FPGA_ID_NOT_CONNECTED ,
[Mea2100Interfaceboard](#) = FPGA_ID_MEA2100_IF ,
[Mea2100Headstage](#) = FPGA_ID_MEA2100_HS ,
[Mea2100STG](#) = FPGA_ID_MEA2100_STG ,
[MultiwellHeadstage](#) = FPGA_ID_HS_MULTIWELL ,
[MultiwellInterfaceboard](#) = FPGA_ID_IF_MULTIWELL ,
[TbsiDacqInterfaceboard](#) = FPGA_ID_TBSI_DACQ_IF ,
[TbsiDacqHeadstage](#) = FPGA_ID_TBSI_DACQ_HS ,

```

CmosMealInterfaceboard = FPGA_ID_CMOS_MEA_IF ,
CmosMeaHeadstage = FPGA_ID_CMOS_MEA_HS ,
Mea2100MultiwellIFB2 = FPGA_ID_MEA2100_MW_IFB2 ,
Me2100Interfaceboard = FPGA_ID_ME2100_IFB ,
Me2100InvivoSignalCollectorUnit = FPGA_ID_ME2100_InvivoSCU ,
Me2100InvitroSignalCollectorUnit = FPGA_ID_ME2100_InvitroSCU ,
Me2100_32XilinxHeadstage = FPGA_ID_ME2100_32_XILINX_HS ,
Me2100_32PICiCE40Headstage = FPGA_ID_ME2100_32_PIC_ICE40_HS ,
Mea2100_256Interfaceboard = FPGA_ID_MEA2100_256_IF ,
Mea2100_256Headstage = FPGA_ID_MEA2100_256_HS ,
W2100Interfaceboard = FPGA_ID_W2100_IF ,
W2100WirelessReceiver = FPGA_ID_W2100_REC ,
W2100WirelessReceiverAnalog = FPGA_ID_W2100_REC_ANALOG ,
Mea2100Mini60PICiCE40Headstage = FPGA_ID_MEA2100MINI60_PIC_ICE40_HS ,
Mea2100BetaScreenHeadstage = FPGA_ID_MEA2100BETASCREEN_HS ,
Me2100UPA32Headstage = FPGA_ID_ME2100UPA32_HS ,
MultiwellMiniHeadstage = FPGA_ID_MULTIWELL_MINI_HS ,
Mea2100Mini120Headstage = FPGA_ID_MEA2100MINI120_HS ,
Mea2100Mini60ECP5Headstage = FPGA_ID_MEA2100MINI60_ECP5_HS ,
eCubeHeadstage = FPGA_ID_ECUBE_HS ,
Me2100Graphene16_32Headstage = FPGA_ID_ME2100_GRAPHENE_16_32_HS ,
GrapheneASICHeadstage = FPGA_ID_GRAPHENE_ASIC_HS ,
WholeCellPatchHeadstage = FPGA_ID_WHOLE_CELL_PATCH_HS ,
InterfaceBoard2 = FPGA_ID_INTERFACEBOARD2 ,
W2100IFB2 = FPGA_ID_W2100_IFB2 ,
CmosmealIFB2 = FPGA_ID_CMOS_MEA_IFB2 ,
Mea2100LiteHeadstage = FPGA_ID_MEA2100_LITE_HS ,
LIH30Interfaceboard = FPGA_ID_LIH30_USB_IF ,
LIH30ADCCtrl = FPGA_ID_LIH30_ADC_CTRL ,
UssingRail = FPGA_ID_USSING_RAIL ,
UssingChamber = FPGA_ID_USSING_CHAMBER ,
IFB2GoldenInterfaceboard = FPGA_ID_IFB2_GOLDEN ,
IFB30GoldenInterfaceboard = FPGA_ID_IFB30_GOLDEN ,
DeviceHasNoHeadstage = FPGA_ID_HAS_NO_HS }
• enum class HeadstageIdEnumNet {
DeviceNotConnected = FPGA_ID_NOT_CONNECTED ,
Mea2100 = FPGA_ID_MEA2100_HS ,
Multiwell = FPGA_ID_HS_MULTIWELL ,
TbsiDacq = FPGA_ID_TBSI_DACQ_HS ,
CmosMea = FPGA_ID_CMOS_MEA_HS ,
InvivoSignalCollectorUnit = FPGA_ID_ME2100_InvivoSCU ,
InvitroSignalCollectorUnit = FPGA_ID_ME2100_InvitroSCU ,
Mea2100_256 = FPGA_ID_MEA2100_256_HS ,
W2100WirelessReceiver = FPGA_ID_W2100_REC ,
W2100WirelessReceiverAnalog = FPGA_ID_W2100_REC_ANALOG ,
Mea2100_Lite = FPGA_ID_MEA2100_LITE_HS ,
LIH30ADCCtrl = FPGA_ID_LIH30_ADC_CTRL ,
DeviceHasNoHeadstage = FPGA_ID_HAS_NO_HS }
• enum class SCU_HeadstageIdEnumNet {
DeviceNotConnected = FPGA_ID_NOT_CONNECTED ,
Me2100_32Xilinx = FPGA_ID_ME2100_32_XILINX_HS ,
Me2100_32PICiCE40 = FPGA_ID_ME2100_32_PIC_ICE40_HS ,
Mea2100Mini60PICiCE40 = FPGA_ID_MEA2100MINI60_PIC_ICE40_HS ,
Mea2100BetaScreen = FPGA_ID_MEA2100BETASCREEN_HS ,
Me2100UPA32 = FPGA_ID_ME2100UPA32_HS ,
MultiwellMini = FPGA_ID_MULTIWELL_MINI_HS ,
Mea2100Mini120 = FPGA_ID_MEA2100MINI120_HS ,
Mea2100Mini60ECP5 = FPGA_ID_MEA2100MINI60_ECP5_HS ,

```

```

eCube = FPGA_ID_ECUBE_HS ,
Me2100Graphene16_32 = FPGA_ID_ME2100_GRAPHENE_16_32_HS ,
GrapheneASIC = FPGA_ID_GRAPHENE_ASIC_HS ,
WholeCellPatch = FPGA_ID_WHOLE_CELL_PATCH_HS ,
DeviceHasNoHeadstage = FPGA_ID_HAS_NO_HS }
• enum class UsbVendorIdEnumNet {
    Unknown = -1 ,
    None = 0 ,
    Renesas = 0x1912 ,
    ASMedia = 0x1b21 ,
    Intel = 0x8086 }
• enum class FilterCalculationDirectionEnumNet {
    DoubleToInt = 0 ,
    IntToDouble = 1 }
• enum class FilterBandEnumNet {
    Unknown = 0 ,
    Lowpass = 1 ,
    Highpass = 2 }
• enum class FilterFamilyEnumNet {
    Unknown = 0 ,
    Bessel = 1 ,
    Butterworth = 2 ,
    RC = 3 }
• enum class FilterTypeEnumNet {
    Hardware = 0 ,
    Software = 1 }
• enum class FilterAttributeEnumNet {
    PreCommaB = 0 ,
    PostCommaB = 1 ,
    CommaPositionB = 2 ,
    PreCommaA = 3 ,
    PostCommaA = 4 ,
    CommaPositionA = 5 }
• enum class AnalogOut_DAC_Range_EnumNet {
    PlusMinus2Comma5Volts = 0 ,
    PlusMinus5Volts = 1 ,
    PlusMinus10Volts = 2 }
• enum class PP_Pump_Mode_Type_EnumNet {
    Manual = 1 ,
    Digital = 2 ,
    Analog = 3 }
• enum class MbcChargingModeEnumNet {
    StorageCharge = 0 ,
    FullCharge = 1 }
• enum class MbcRatedCapacityEnumNet {
    rc30mAh = 0 ,
    rc100mAh = 1 ,
    rc200mAh = 2 ,
    rc300mAh = 3 ,
    rcGreater300mAh = 4 }
• enum class MbcChannelStateEnumNet {
    csIdleNoBattery = 0 ,
    csIdleChargeFinished = 1 ,
    csCapacityTestPrecharge = 2 ,
    csCapacityTestDischarge = 3 ,
    csRefreshBattery = 4 ,
    csCharge = 5 ,

```

```

    csDischarge = 6 ,
    csError = 7 }
• enum class PulseGenerator_Mode_EnumNet {
    Off = 0 ,
    AlwaysOn = 1 ,
    Gated_Low_Active = 2 ,
    Gated_High_Active = 3 }
• enum class LIH30_ADC_Channel_EnumNet {
    User_ADC_0 = 0 ,
    User_ADC_1 = 1 ,
    User_ADC_2 = 2 ,
    User_ADC_3 = 3 ,
    User_ADC_4 = 4 ,
    Test_ADC_EPC10 = 5 ,
    ModulA_ADC0 = 6 ,
    ModulA_ADC1 = 7 ,
    ModulA_ADC2 = 8 ,
    ModulA_ADC3 = 9 ,
    ModulB_ADC0 = 10 ,
    ModulB_ADC1 = 11 ,
    ModulB_ADC2 = 12 ,
    ModulB_ADC3 = 13 ,
    ModulC_ADC0 = 14 ,
    ModulC_ADC1 = 15 ,
    ModulC_ADC2 = 16 ,
    ModulC_ADC3 = 17 ,
    ModulD_ADC0 = 18 ,
    ModulD_ADC1 = 19 ,
    ModulD_ADC2 = 20 ,
    ModulD_ADC3 = 21 }
• enum class LIH30_DAC_Channel_EnumNet {
    User_DAC_0 = 0 ,
    User_DAC_1 = 1 ,
    User_DAC_2 = 2 ,
    Test_DAC_EPC10 = 3 ,
    ModulA_DAC0 = 4 ,
    ModulA_DAC1 = 5 ,
    ModulB_DAC0 = 6 ,
    ModulB_DAC1 = 7 ,
    ModulC_DAC0 = 8 ,
    ModulC_DAC1 = 9 ,
    ModulD_DAC0 = 10 ,
    ModulD_DAC1 = 11 }
• enum class LIH30_EPC10_Bus_EnumNet {
    A = 0 ,
    B = 1 }
• enum class W2100_Accel_Gyro_Select_EnumNet {
    Off = 0 ,
    GyroOnly = 1 ,
    AccelOnly = 2 ,
    Both = 3 }
    enumerates the accelerometer configuration on the W2100 device
• enum class WvcValveModeEnumNet {
    Manual = WVC_VALVE_MODE_MANUAL ,
    Digital = WVC_VALVE_MODE_DIGITAL ,
    Analog = WVC_VALVE_MODE_ANALOG ,
    Table = WVC_VALVE_MODE_TABLE }
    enumerates Wvc valve mode

```

- enum class [WvcDisplayModeEnumNet](#) {
[Work](#) = WVC_DISPLAY_MODE_WORK ,
[PC](#) = WVC_DISPLAY_MODE_PC ,
[Settings](#) = WVC_DISPLAY_MODE_SETTINGS ,
[TouchTest](#) = WVC_DISPLAY_MODE_TOUCH_TEST }
enumerates Wvc display mode
- enum class [PortDirectionEnumNet](#) {
[Output](#) = 0 ,
[Input](#) = 1 }
enumerates a port direction
- enum class [StimulationLayoutConfigurationEnumNet](#) {
[SingleWell](#) = 1 ,
[SixWell](#) = 2 ,
[NineWell](#) = 3 }
enumerates the layout configuration for the MEA2100-256 device
- enum class [ReferenceElectrodeSwitchPositionEnumNet](#) {
[off](#) = 0 ,
[Ref8](#) = 1 ,
[Ref16](#) = 2 ,
[Ref24](#) = 3 ,
[Ref32](#) = 4 }
enumerates the possible positions of the reference electrode switch of the ME2100 device
- enum class [ReferenceElectrodeModeEnumNet](#) {
[SubtractionOff](#) = 0 ,
[SubtractFromAllOther](#) = 1 ,
[SubtractFromReferenceElectrodeOnly](#) = 2 ,
[SubtractFromAll](#) = 3 }
enumerates the electrode subtraction modes
- enum class [DigitalDatastreamEnableEnumNet](#) {
[None](#) = 0x0000 ,
[Mux](#) = 0x0001 ,
[MuxOtherDevice](#) = 0x0002 ,
[DigitalInReserverd](#) = 0x0004 ,
[DigitalIn](#) = 0x0008 ,
[DigitalOut](#) = 0x0010 ,
[DigitalOutReserved](#) = 0x0020 ,
[RegisterLow](#) = 0x0040 ,
[RegisterHigh](#) = 0x0080 ,
[FeedbackLow](#) = 0x0100 ,
[FeedbackHigh](#) = 0x0200 ,
[Aux](#) = 0x0400 ,
[PeriodicPulse](#) = 0x0800 ,
[DigOutStim](#) = 0x1000 ,
[Hs1Digital](#) = 0x00008000 ,
[Hs1Trigger](#) = 0x00010000 ,
[Hs1SidebandLow](#) = 0x00020000 ,
[Hs1SidebandHigh](#) = 0x00040000 ,
[Hs2Digital](#) = 0x00800000 ,
[Hs2Trigger](#) = 0x01000000 ,
[Hs2SidebandLow](#) = 0x02000000 ,
[Hs2SidebandHigh](#) = 0x04000000 }
enumerates the streams available as digital datastream
- enum class [IoVoltageEnumNet](#) {
[Voltage_3V3](#) = IFB2_IO_VOLTAGE_3V3 ,
[Voltage_5V0](#) = IFB2_IO_VOLTAGE_5V0 }
enumerates the I/O Voltages available on the IFB2

- enum class `EnSTG200x_STATUS` {
`OK` ,
`NOT_CONNECTED` ,
`DEVICE_NOT_FOUND` }

Functions

- public delegate void `OnMcsUsbDeviceState` (`usbSetupPacket_t`[^] request)
- private delegate void `OnMcsUsbDeviceStateCallback` (IntPtr pThis, uint32_t size, IntPtr buffer)
- public delegate void `OnUpdateFirmwareStatusChange` (String[^])
- public delegate void `OnUpdateFirmwareProgress` (int)
- public delegate void `OnDeviceArrivalRemoval` (`CMcsUsbListEntryNet`[^] entry)
Delegate to show a device arrival or removal.
- public delegate void `OnStgPollStatus` (unsigned int status, `StgStatusNet`[^] stgStatusNet, array< int >[^] index_list)
- public delegate void `OnMwPollStatus` (unsigned int CurrentTemp, unsigned int PlateState, unsigned int SwitchState)
- public delegate void `RoboStatusEventDelegate` (array< unsigned char >[^] buffer)
- public delegate void `OnStg200xDataHandler` (uint32_t trigger)
- public delegate void `OnStg200xErrorHandler` ()
- public delegate void `OnChannelData` (`CMcsUsbDacqNet`[^] dacq, int CbHandle, int numFrames)
- public delegate void `OnError` (String[^] msg, int action)

10.2.1 Enumeration Type Documentation

10.2.1.1 AdapterTypeEnumNet enum AdapterTypeEnumNet [strong]

Enumerates the adapter type of the MEA2100 device.

Enumerator

None	
MEA60	
MEA2x60	
MEA120	
MEA32	
MEA2x32	
Multiwell96	
WirelessTestAdapter	
MEA252	
MEA_2_252_2	
MEA_2_252_2_6Well	
MEA_2_252_2_9Well	
MEA_2_252_2_Test	
TBSI_5	
TBSI_15	
TBSI_31	
TBSI_63	
TBSI_127	
TBSI_Reserved	
Ci4600Intan	
Unknown	
NotApplicable	

10.2.1.2 AnalogOut_DAC_Range_EnumNet enum [AnalogOut_DAC_Range_EnumNet](#) [strong]

Enumerator

PlusMinus2Comma5Volts	
PlusMinus5Volts	
PlusMinus10Volts	

10.2.1.3 AnalogSourceEnumNet enum [AnalogSourceEnumNet](#) [strong]

Enumerates the analog source of the MEA2100 device.

Enumerator

AnalogSource_HS1	
AnalogSource_HS2	
AnalogSource_IF	

10.2.1.4 AnalogUnitEnumNet enum [AnalogUnitEnumNet](#) [strong]

Enumerator

Unknown	
Volt	
Ampere	
Kelvin	

10.2.1.5 CFirmwareDestinationNet enum [CFirmwareDestinationNet](#) [strong]

Enumerates the destination processor for the firmware.

Enumerator

FPGA_NORMAL	
DSP	The DSP.
USB	The USB controller.
MCU1	The DSP on the MEA2100 system.
MCSBUS1	
MCSBUS2	

Enumerator

MCSBUS3	
MCSBUS4	
MCSBUS5	
MCSBUS6	
MCSBUS7	
MCSBUS8	
MCSBUS9	
MCSBUS10	
MCSBUS11	
MCSBUS12	
MCSBUS13	
MCSBUS14	
MCSBUS15	
MCSBUS0	
BUSNUMBER1	
BUS1MCSBUS1	
BUS1MCSBUS2	
BUSNUMBER2	
BUS2MCSBUS1	
BUS2MCSBUS2	
PIC	
PIC2	
PIC3	
PIC4	
PIC5	
PIC6	
PIC7	
PIC8	
PIC9	
PIC10	
PIC11	
PIC12	
ChannelPIC	
Bootstrap	
BootstrapOtherCypress	
ALTERA	
FPGA2	
FPGA3	
FPGA4	
FPGA5	
FPGA6	
FPGA7	
FPGA8	
FPGA9	
FPGA10	
FPGA11	
FPGA12	
FPGA13	
FPGA14	
FPGA15	

Enumerator

FPGA16	
FPGA_GOLD	
ALTERA_GOLD	
FPGA2_GOLD	
FPGA3_GOLD	
FPGA4_GOLD	
FPGA5_GOLD	
FPGA6_GOLD	
FPGA7_GOLD	
FPGA8_GOLD	
FPGA9_GOLD	
FPGA10_GOLD	
FPGA11_GOLD	
FPGA12_GOLD	
FPGA13_GOLD	
FPGA14_GOLD	
FPGA15_GOLD	
FPGA16_GOLD	
FPGA_BASE	
ALTERA_BASE	
FPGA2_BASE	
FPGA3_BASE	
FPGA4_BASE	
FPGA5_BASE	
FPGA6_BASE	
FPGA7_BASE	
FPGA8_BASE	
FPGA9_BASE	
FPGA10_BASE	
FPGA11_BASE	
FPGA12_BASE	
FPGA13_BASE	
FPGA14_BASE	
FPGA15_BASE	
FPGA16_BASE	
FPGA_BOOTSTRAP	
ALTERA_BOOTSTRAP	
DEST_TARGET1	
DEST_TARGET2	
DEST_TARGET3	
DEST_TARGET4	
DEST_TARGET5	
DEST_TARGET6	
DEST_TARGET7	
DEST_TARGET8	
DEST_TARGET9	
DEST_TARGET10	
DEST_TARGET11	
DEST_TARGET12	

Enumerator

DEST_TARGET13	
DEST_TARGET14	
DEST_TARGET15	
DEST_TARGET_MASK	
DEST_FX3_TARGET_MASK	
ALTERA_TARGET1	
ALTERA_TARGET2	
ALTERA_TARGET3	
USB_TARGET1	
USB_TARGET2	
USB_TARGET3	
UnknownDest	

10.2.1.6 CMOSMeaBathModeEnumNet `enum CMOSMeaBathModeEnumNet [strong]`

Enumerations of CMOS MEA Bath Mode

Enumerator

Ground	
Stimulation	
CurrentMeasure	

10.2.1.7 CMOSMeaHeadstage1NCBathCurrentEnumNet `enum CMOSMeaHeadstage1NCBathCurrentEnumNet [strong]`

Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

Enumerator

NCBathCurrent	
---------------	--

10.2.1.8 CMOSMeaHeadstage1NCCol2CurrentEnumNet `enum CMOSMeaHeadstage1NCCol2CurrentEnumNet [strong]`

Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

Enumerator

NCCol2Current	
---------------	--

10.2.1.9 CMOSMeaHeadstage1NChipTempEnumNet enum [CMOSMeaHeadstage1NChipTempEnumNet](#) [strong]

Enumerations of CMOS MEA HS Temperature Monitoring Channels Group Bitmask

Enumerator

NChipTemperature	
------------------	--

10.2.1.10 CMOSMeaHS1SidebandEnumNet enum [CMOSMeaHS1SidebandEnumNet](#) [strong]

Enumerations of CMOS MEA HS STG Sideband Channels Group Bitmask

Enumerator

SBSVector1	
SBSVector2	
SBSVector3	
SBSVector4	

10.2.1.11 CMOSMeaHS1TriggerStatusEnumNet enum [CMOSMeaHS1TriggerStatusEnumNet](#) [strong]

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

Enumerator

TriggerStatus1	
TriggerStatus2	
TriggerStatus3	
TriggerStatus4	

10.2.1.12 CMOSMeaIFDigChannelEnumNet enum [CMOSMeaIFDigChannelEnumNet](#) [strong]

Enumerations of CMOS MEA IF Digital Channels Group Bitmask

Enumerator

DigitalMux	
DigitalInPort	
DigitalOutReg	

Enumerator

FeedbackReg	
DigitalReg	
AuxPort	

10.2.1.13 CMOSMeaInterfaceADCEnumNet enum [CMOSMeaInterfaceADCEnumNet](#) [strong]

Enumerations of CMOS MEA IF Analog Channels Group Bitmask

Enumerator

IFChannel1	
IFChannel2	
IFChannel3	
IFChannel4	
IFChannel5	
IFChannel6	
IFChannel7	
IFChannel8	

10.2.1.14 CMOSMeaPacketFrameContextGroupEnumNet enum [CMOSMeaPacketFrameContextGroupEnumNet](#) [strong]

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

Enumerator

SOFAndCTRLword	
ChecksumAndPacketCounter	
Timestamp	
EOFAndCRC	

10.2.1.15 CMOSMeaSTG1DACSignalEnumNet enum [CMOSMeaSTG1DACSignalEnumNet](#) [strong]

Enumerations of CMOS MEA DAC Stimulation Channels Group Bitmask

Enumerator

DAC1Channel	
DAC2Channel	
DAC3Channel	
DAC4Channel	

10.2.1.16 CMOSMeaValueUnitEnumNet enum [CMOSMeaValueUnitEnumNet](#) [strong]

Enumerations of CMOS MEA Units of Values in Data stream

Enumerator

NoUnit	
NanoVolt	
PicoAmpere	
NanoAmpere	
MicroAmpere	
MilliDegreeCelsius	

10.2.1.17 DacqGroupChannelEnumNet enum [DacqGroupChannelEnumNet](#) [strong]

Enumerates the Channel Groups of Datastream

Enumerator

HeadstageElectrodeGroup	
InterfaceADCGroup	
DSPDataGroup	
Headstage1NCBathCurrentGroup	
Headstage1NCCol2CurrentGroup	
Headstage1NChipTempGroup	
STG1DACSignalGroup	
LIH30UserADCGroup	
LIH30TestADCGroup	
LIH30ADCModulesGroup	
IFDigChannelsGroup	
STG1SidebandsGroup	
STG1TriggerStatusGroup	
DACQ1DigitalGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.18 DacqMeaGroupTypeEnumNet enum [DacqMeaGroupTypeEnumNet](#) [strong]

Enumerations of CMOS MEA Groups to detect wether it is an Analog, Digital or Frame Context Group

Enumerator

AnalogGroup	
DigitalGroup	
FrameContextGroup	

10.2.1.19 DataModeEnumNet enum [DataModeEnumNet](#) [strong]

Enumerates the data mode of the device, either 16, 24 or 32 bit, can be signed or unsigned.

Enumerator

Unsigned_16bit	
Unsigned_24bit	
Unsigned_32bit	
Signed_16bit	
Signed_24bit	
Signed_32bit	

10.2.1.20 DeviceEnumNet enum [DeviceEnumNet](#) [strong]

Enumerates the group of MCS devices to connect to.

Enumerator

MCS_DEVICE_ANY	To connect to any MCS device.
MCS_GENERIC_DEVELOPMENT_DEVICE	Please use this only for MCS internal development.
MCS_DEVICE_USB	To connect to any MCS USB device.
MCS_MCCARD_DEVICE	Connect to an MC_Card.
MCS_STG_DEVICE	Connect to an MCS device with STG capability.
MCS_MC_STIMULUS_DEVICE	Devices which should be accessible from MC_Stimulus.
MCS_MEASUSB_DEVICE	Connect to an MCS MeaUsb device.
MCS_MEA_DEVICE	Connect to an MCS MeaUsb device.
MCS_OCTOPOT_DEVICE	Connect to an MCS Octopot device.
MCS_TERSENS_DEVICE	Connect to an MCS Tersens device.
MCS_PGA_DEVICE	Connect to an MCS PGA device.
MCS_PCX_DEVICE	Connect to an MCS PCX device.
MCS_TCX_DEVICE	Connect to an MCS TCX device.
MCS_FCX_DEVICE	Connect to an MCS FCX device.
MCS_RETINA_LED_DEVICE	Connect to an MCS RetineLed device.
MCS_MEA_SWITCH_DEVICE	Connect to an MCS Mea Switch device.
MCS_MEA_IMPEDANCE_DEVICE	Connect to an MCS Mea Impedance device.
MCS_CHANNELTEST_DEVICE	Connect to an MCS ChannelTest device.
MCS_SW2TO64_DEVICE	Connect to an MCS SW2TO64 device.

Enumerator

MCS_RETINA_AMS_DONGLE	Connect to an MCS Retina AMS Dongle (Radio device)
MCS_PATHIDENT_DEVICE	Connect to an MCS Pathident device.
MCS_ROBO_DEVICE	Connect to an MCS Robo Platform device.
MCS_ROBOOCYTE2_DEVICE	Connect to an MCS Roboocyte2 device.
MCS_ROBOINJECT_DEVICE	Connect to an MCS RoboInject device.
MCS_HICLAMP_DEVICE	Connect to an MCS HiClamp device.
MCS_PATCHSERVER_DEVICE	Connect to an MCS PatchServer device.
MCS_ENCAPSULATOR_DEVICE	Connect to an MCS Encapsulator device.
MCS_MEASURETABLE_DEVICE	Connect to an MCS Encapsulator device.
MCS_FYI_DEVICE	Connect to an MCS FYI device.
MCS_HLA_DEVICE	Connect to an MCS HLA device.
MCS_PPS_DEVICE	Connect to an MCS PPS device.
MCS_PPS5_DEVICE	Connect to an MCS PPS5 device.
MCS_OKUVISION_STIMULATOR_DEVICE	Connect to an Okuvision Stimulator device.
MCS_NF_GEN_DEVICE	Connect to an MCS NF-Gen device.
MCS_SAFEIS_DEVICE	Connect to an MCS SafeIS device.
MCS_PERISTALTIC_PUMP_DEVICE	Connect to an MCS PeristalticPump device.
MCS_EXTERN_BC_TESTER_DEVICE	Connect to an ExternBCTester device.
MCS_EXTERN_D_TESTER_DEVICE	Connect to an ExternDTester device.
MCS_SOFTWARE_DONGLE_DEVICE	Connect to an Software Dongle device.
MCS_MEA_CLEAN_DEVICE	Connect to a MEA Clean device.
MCS_MEA_COAT_DEVICE	Connect to a MEA Clean device.
MCS_SMARTIMPLANT_DEVICE	Connect to a SmartImplant device.
MCS_MBC08_DEVICE	Connect to a MultiBatteryCharger device.
MCS_PEDOTER_DEVICE	Connect to a Pedoter device.
MCS_PPC_DEVICE	Connect to a PPC device.
WARNER_VALVE_CONTROL_DEVICE	Connect to a Warner Valve Control device.
WARNER_USSING_DEVICE	Connect to a Warner Valve Control device.
HEKA_LIH3_DEVICE	Connect to a HEKA LIH3 device.
ALA_VC3_DEVICE	Connect to an ALA VC3 Valve Commander.
MCS_DEVICE_USB_CYPRESS	Connect to a Cypress USB device.

10.2.1.21 DigitalDatastreamEnableEnumNet `enum DigitalDatastreamEnableEnumNet` [strong]

enumerates the streams available as digital datastream

Enumerator

None	No digital datastream.
Mux	16 bits from the standard MUX datastream.
MuxOtherDevice	The 16 bits of the standard MUX datastream used by the other virtual device.
DigitalInReserverd	The lower 16 bits of the Digital IN port, these ports are on the device by default used as Digital OUT, thus not available as Digital IN.
DigitalIn	The upper 16 bits of the Digital IN port, use when Digital IN datastream is needed.
DigitalOut	The lower 16 bits of the Digital OUT port, use when Digital OUT datastream is needed.

Enumerator

DigitalOutReserved	The upper 16 bits of the Digital OUT port, these ports are on the device by default used as Digital IN, thus not available as Digital OUT.
RegisterLow	The lower 16 bits of the digital register.
RegisterHigh	The upper 16 bits of the digital register.
FeedbackLow	The lower 16 bits of the realtime feedback register.
FeedbackHigh	The upper 16 bits of the realtime feedback register.
Aux	The 2 bits of the AUX port.
PeriodicPulse	The 8 bits of the Periodic Pulse Generator (Video-Sync).
DigOutStim	The 16 bits of the Digital Out Stimulator.
Hs1Digital	Headstage 1 digital signals.
Hs1Trigger	Headstage 1 trigger signals.
Hs1SidebandLow	Headstage 1 lower 16 bits of sideband data.
Hs1SidebandHigh	Headstage 1 upper 16 bits of sideband data.
Hs2Digital	Headstage 2 digital signals.
Hs2Trigger	Headstage 2 trigger signals.
Hs2SidebandLow	Headstage 2 lower 16 bits of sideband data.
Hs2SidebandHigh	Headstage 2 upper 16 bits of sideband data.

10.2.1.22 DigitalSourceEnumNet `enum DigitalSourceEnumNet [strong]`

Enumerates the digital source of the MEA2100 device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Trigger3Status	
HS1Trigger4Status	
HS1Trigger5Status	
HS1Trigger6Status	
HS1Sideband1	
HS1Sideband2	
HS1Sideband3	
HS1Sideband4	
HS1Sideband5	
HS1Sideband6	
HS2Trigger1Status	
HS2Trigger2Status	

Enumerator

HS2Trigger3Status	
HS2Trigger4Status	
HS2Trigger5Status	
HS2Trigger6Status	
HS2Sideband1	
HS2Sideband2	
HS2Sideband3	
HS2Sideband4	
HS2Sideband5	
HS2Sideband6	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
LastPosition	

10.2.1.23 DigitalStimulatorTriggerEventEnumNet enum [DigitalStimulatorTriggerEventEnumNet](#) [strong]

Enumerates start/stop event for DigOut/DigStim trigger. /summary>

Enumerator

Start	
Stop	

10.2.1.24 DigitalStimulatorTriggerSlopeEnumNet enum [DigitalStimulatorTriggerSlopeEnumNet](#) [strong]

Enumerates start/stop conditions for DigOut/DigStim trigger. /summary>

Enumerator

Falling	
Rising	

10.2.1.25 DigitalTargetEnumNet enum [DigitalTargetEnumNet](#) [strong]

Enumerates the Digital Targets for Digital Sources

Enumerator

Digout	
Digstream	
DacqTrigger	
StgTrigger	
StgListModeTrigger	
DigOutStimulatorStartTrigger	
DigOutStimulatorStopTrigger	
DigStreamToReceiver	

10.2.1.26 ElectrodeDacMuxEnumNet enum [ElectrodeDacMuxEnumNet](#) [strong]

Enumerates the setting of the Stimulation DAC Multiplexer.

Enumerator

Ground	Connect the electrode to Ground while stimulation is active.
Stg1	Connect the electrode tp STG 1 while stimulation is active.
Stg2	Connect the electrode tp STG 2 while stimulation is active.
Stg3	Connect the electrode tp STG 3 while stimulation is active.

10.2.1.27 ElectrodeModeEnumNet enum [ElectrodeModeEnumNet](#) [strong]

Enumerates the mode of each electrode, can be automatic or manual. In automatic mode, the blanking of the electrode is controlled by the sideband signal, in manual mode, the stimulation configuration is independant of the sideband signal.

Enumerator

emAutomatic	
emManual	

10.2.1.28 enCMosMeaChipType enum [enCMosMeaChipType](#) [strong]

Enumerator

unknown	
nMos16LV	
nMos32LV	
nMos36LN	
nMos64LN	

10.2.1.29 EnSTG200x_STATUS enum [EnSTG200x_STATUS](#) [strong]

Enumerator

OK	
NOT_CONNECTED	
DEVICE_NOT_FOUND	

10.2.1.30 FilterAttributeEnumNet enum [FilterAttributeEnumNet](#) [strong]

Enumerator

PreCommaB	
PostCommaB	
CommaPositionB	
PreCommaA	
PostCommaA	
CommaPositionA	

10.2.1.31 FilterBandEnumNet enum [FilterBandEnumNet](#) [strong]

Enumerator

Unknown	
Lowpass	
Highpass	

10.2.1.32 FilterCalculationDirectionEnumNet enum [FilterCalculationDirectionEnumNet](#) [strong]

Enumerator

DoubleToInt	
IntToDouble	

10.2.1.33 FilterFamilyEnumNet enum [FilterFamilyEnumNet](#) [strong]

Enumerator

Unknown	
Bessel	
Butterworth	
RC	

10.2.1.34 FilterTypeEnumNet enum `FilterTypeEnumNet` [strong]

Enumerator

Hardware	
Software	

10.2.1.35 FpgaldEnumNet enum `FpgaIdEnumNet` [strong]

Enumerator

DeviceNotConnected	
Mea2100Interfaceboard	
Mea2100Headstage	
Mea2100STG	
MultiwellHeadstage	
MultiwellInterfaceboard	
TbsiDacqInterfaceboard	
TbsiDacqHeadstage	
CmosMeaInterfaceboard	
CmosMeaHeadstage	
Mea2100MultiwellIFB2	
Me2100Interfaceboard	
Me2100InvivoSignalCollectorUnit	
Me2100InvitroSignalCollectorUnit	
Me2100_32XilinxHeadstage	
Me2100_32PICI40Headstage	
Mea2100_256Interfaceboard	
Mea2100_256Headstage	
W2100Interfaceboard	
W2100WirelessReceiver	
W2100WirelessReceiverAnalog	
Mea2100Mini60PICI40Headstage	
Mea2100BetaScreenHeadstage	
Me2100UPA32Headstage	
MultiwellMiniHeadstage	
Mea2100Mini120Headstage	

Enumerator

Mea2100Mini60ECP5Headstage	
eCubeHeadstage	
Me2100Graphene16_32Headstage	
GrapheneASICHeadstage	
WholeCellPatchHeadstage	
InterfaceBoard2	
W2100IFB2	
CmosmealFB2	
Mea2100LiteHeadstage	
LIH30Interfaceboard	
LIH30ADCCtrl	
UssingRail	
UssingChamber	
IFB2GoldenInterfaceboard	
IFB30GoldenInterfaceboard	
DeviceHasNoHeadstage	

10.2.1.36 HeadstageIdEnumNet enum [HeadstageIdEnumNet](#) [strong]

Enumerator

DeviceNotConnected	
Mea2100	
Multiwell	
TbsiDacq	
CmosMea	
InvivoSignalCollectorUnit	
InvitroSignalCollectorUnit	
Mea2100_256	
W2100WirelessReceiver	
W2100WirelessReceiverAnalog	
Mea2100_Lite	
LIH30ADCCtrl	
DeviceHasNoHeadstage	

10.2.1.37 IoVoltageEnumNet enum [IoVoltageEnumNet](#) [strong]

enumerates the I/O Voltages available on the IFB2

Enumerator

Voltage_3V3	
Voltage_5V0	

10.2.1.38 LIH30_ADC_Channel_EnumNet enum [LIH30_ADC_Channel_EnumNet](#) [strong]

Enumerator

User_ADC_0	
User_ADC_1	
User_ADC_2	
User_ADC_3	
User_ADC_4	
Test_ADC_EPC10	
ModulA_ADC0	
ModulA_ADC1	
ModulA_ADC2	
ModulA_ADC3	
ModulB_ADC0	
ModulB_ADC1	
ModulB_ADC2	
ModulB_ADC3	
ModulC_ADC0	
ModulC_ADC1	
ModulC_ADC2	
ModulC_ADC3	
ModulD_ADC0	
ModulD_ADC1	
ModulD_ADC2	
ModulD_ADC3	

10.2.1.39 LIH30_DAC_Channel_EnumNet enum [LIH30_DAC_Channel_EnumNet](#) [strong]

Enumerator

User_DAC_0	
User_DAC_1	
User_DAC_2	
Test_DAC_EPC10	
ModulA_DAC0	
ModulA_DAC1	
ModulB_DAC0	
ModulB_DAC1	
ModulC_DAC0	
ModulC_DAC1	
ModulD_DAC0	
ModulD_DAC1	

10.2.1.40 LIH30_EPC10_Bus_EnumNet enum [LIH30_EPC10_Bus_EnumNet](#) [strong]

Enumerator

A	
B	

10.2.1.41 MbcChannelStateEnumNet enum [MbcChannelStateEnumNet](#) [strong]

Enumerator

csIdleNoBattery	
csIdleChargeFinished	
csCapacityTestPrecharge	
csCapacityTestDischarge	
csRefreshBattery	
csCharge	
csDischarge	
csError	

10.2.1.42 MbcChargingModeEnumNet enum [MbcChargingModeEnumNet](#) [strong]

Enumerator

StorageCharge	
FullCharge	

10.2.1.43 MbcRatedCapacityEnumNet enum [MbcRatedCapacityEnumNet](#) [strong]

Enumerator

rc30mAh	
rc100mAh	
rc200mAh	
rc300mAh	
rcGreater300mAh	

10.2.1.44 McsBusTypeEnumNet enum [McsBusTypeEnumNet](#) [strong]

Enumerates the bus to use, either USB, PCI or any

Enumerator

MCS_ANY_BUS	
MCS_UNDEFINED_BUS	
MCS_USB_BUS	
MCS_PCI_BUS	

10.2.1.45 McsUsbSpeedEnumNet enum [McsUsbSpeedEnumNet](#) [strong]

Enumerates the current connection speed of the device

Enumerator

LowSpeed	
FullSpeed	
HighSpeed	
SuperSpeed	
UnknownSpeed	

10.2.1.46 MEA2100_256DacqGroupChannelEnumNet enum [MEA2100_256DacqGroupChannelEnumNet](#) [strong]

Enumerates the MEA2100-256 Channel Groups of Datastream

Enumerator

HS1ElectrodeGroup	
HS2ElectrodeGroup	
InterfaceADCGroup	
STG1DACSignalGroup	
STG2DACSignalGroup	
DSPAnalogGroup	
DSPDigitalGroup	
IFDigChannelsGroup	
STG1TriggerStatusGroup	
STG1SidebandsGroup	
STG2TriggerStatusGroup	
STG2SidebandsGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.47 MEA2100_256DigitalSourceEnumNet enum [MEA2100_256DigitalSourceEnumNet](#) [strong]

Enumerates the digital source of the MEA2100-256 device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
DeviceRunStatus	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Trigger3Status	
HS1Trigger4Status	
HS1Trigger5Status	
HS1Trigger6Status	
HS1Trigger7Status	
HS1Trigger8Status	
HS1Trigger9Status	
HS1Trigger10Status	
HS1Trigger11Status	
HS1Trigger12Status	
HS1Trigger13Status	
HS1Trigger14Status	
HS1Trigger15Status	
HS1Trigger16Status	
HS1Trigger17Status	
HS1Trigger18Status	
HS1Sideband1	
HS1Sideband2	
HS1Sideband3	
HS1Sideband4	
HS1Sideband5	
HS1Sideband6	
HS1Sideband7	
HS1Sideband8	
HS1Sideband9	
HS1Sideband10	
HS1Sideband11	
HS1Sideband12	
HS1Sideband13	
HS1Sideband14	
HS1Sideband15	
HS1Sideband16	
HS1Sideband17	
HS1Sideband18	

Enumerator

HS2Trigger1Status	
HS2Trigger2Status	
HS2Trigger3Status	
HS2Trigger4Status	
HS2Trigger5Status	
HS2Trigger6Status	
HS2Trigger7Status	
HS2Trigger8Status	
HS2Trigger9Status	
HS2Trigger10Status	
HS2Trigger11Status	
HS2Trigger12Status	
HS2Trigger13Status	
HS2Trigger14Status	
HS2Trigger15Status	
HS2Trigger16Status	
HS2Trigger17Status	
HS2Trigger18Status	
HS2Sideband1	
HS2Sideband2	
HS2Sideband3	
HS2Sideband4	
HS2Sideband5	
HS2Sideband6	
HS2Sideband7	
HS2Sideband8	
HS2Sideband9	
HS2Sideband10	
HS2Sideband11	
HS2Sideband12	
HS2Sideband13	
HS2Sideband14	
HS2Sideband15	
HS2Sideband16	
HS2Sideband17	
HS2Sideband18	
LastPosition	

10.2.1.48 MeaLayoutEnumNet enum [MeaLayoutEnumNet](#) [strong]

Enumerates the MEA layout of the MEA2100 device.

Enumerator

mlUnknown	
mlMEA60	

10.2.1.49 MultiwellPlateTypeEnumNet enum `MultiwellPlateTypeEnumNet` [strong]

Enumerator

Plate_Dummy	
Plate_24W700_100FMA	
Plate_24W030MGA	
Plate_72W500_100PMA	
Plate_72W500_100FMA	
Plate_24W700_100FMB	
Plate_96W700_100FMA	
Plate_96W300_80_1152FMA	
Plate_96W400_80_1152FMB	
Plate_24W300_30_1152GBA	
Plate_24W700_100FMC	
Plate_96W700_100FMB	
Plate_96W700_100GBC	
Plate_96W700_100GBD	
Plate_24W700_100PBA	
Plate_Dummy_126	
Plate_24W300_30GMA	
Plate_96W700_100GMA	
Plate_24W300_30GBA	
Plate_96W700_100GBA	
Plate_24W300_30GBB	
Plate_96W700_100GBB	
No_Plate	

10.2.1.50 PatchServAdcModeEnumNet enum `PatchServAdcModeEnumNet` [strong]

Enumerator

Normal	
CatchAmp	

10.2.1.51 PlateClampEnumNet enum `PlateClampEnumNet` [strong]

Enumerator

Close	
Open	
Stop	

10.2.1.52 PlateClampLockEnumNet enum [PlateClampLockEnumNet](#) [strong]

Enumerator

Lock	
Unlock	

10.2.1.53 PortDirectionEnumNet enum [PortDirectionEnumNet](#) [strong]

enumerates a port direction

Enumerator

Output	
Input	

10.2.1.54 PP_Pump_Mode_Type_EnumNet enum [PP_Pump_Mode_Type_EnumNet](#) [strong]

Enumerator

Manual	
Digital	
Analog	

10.2.1.55 ProductIdEnumNet enum [ProductIdEnumNet](#) [strong]

Enumerates the group of MCS devices to connect to.

Enumerator

Any	
None	
ALA_VC3	
Cypress_FX1	
Cypress_FX2	
Cypress_FX3	
MC_Card	
Campden_Ci4600EphysVideoDataIntegrator	
HekaLIH30	

Enumerator

HekaEPC10Single	
HekaEPC10Double	
HekaEPC10Triple	
HekaEPC10Quadro	
HekaLIH406	
HekaLIH816	
HekaITEV100	
HekaPG610	
HekaPG611	
HekaPG612	
HekaPG618	
HekaPG690	
HekaEPCLite	
STG	
Octopot	
Tersens	
Dotriapot	
HLA	
STG400x	
STG4002	
STG4004	
STG4008	
STG400x_opto	
STG4002_opto	
STG4004_opto	
STG4008_opto	
STG5	
STG3008_FA	
MultiwellOptoStim	
Generic	
PGA	
PCX	
TCX	
FCX	
FCB	
TC01	
TC02	
Retina_LED	
AMS_Dongle	
Okuvision_Stimulator	
ExternBCTester	
Triggerbox_IMS	
Triggerbox_AMS	
Triggerbox_AMS3	
ExternDTester	
FunkDongleS	
ExternSTester	
DongleS	
Triggerbox_R5	

Enumerator

MEA_Switch	
MEA_Impedance	
ChannelTest	
Sw2to64	
PeristalticPump	
MEA_Switch_2_1	
MEA_Switch_4_2	
PPS4plus1	
PPS5	
PPS2	
PPS5_DIG	
MEA_Clean	
MEA_Coat	
Multiwell_ICC	
MBC08	
PPC	
MEA1060	
MEA_Sanofi	
ME256	
ME128	
ME64	
ME32	
ME16	
MEA2100_Mini_Usb_develop	
MEA256	
MEA2100	
MEA2100_32	
MEA2100_Lite	
Multiwell	
MEA2100_256	
ME2100	
MEA2100BetaScreen	
MEA2100_Mini	
TBSI_Dacq	
Multiwell_MEA_Mini	
Whole_Cell_Patch	
eCube	
Graphene_ASIC	
GE2100	
Multiboot	
WPA8	
WPA4	
WPA16	
WPA32	
W2100	
NeuroChip	
UsbTest	
SoftwareDongle	
PathIdent	
NF_Gen	

Enumerator

SafeIS	
Encapsulator	
NeurochipConfig	
MeasureTable	
Roboocyte2	
RoboInject	
HiClamp	
PatchServer	
Dilutor	
HiClamp4Uart	
IM16S16KRA	
IM64KRB	
IS32KRA	
IM64KRC	
IM16S8KRA	
IM16KRC	
SmartImplant	
PositionImp	
PositionBase	
PositionIICentralUnit	
PositionIIBase	
GrapheneProjectTestDevice	
Pos900	
Neptun	
Warner_Valve_Control	
Warner_TEER_Machine	
Warner_Ussing	

10.2.1.56 PulseGenerator_Mode_EnumNet enum [PulseGenerator_Mode_EnumNet](#) [strong]

Enumerator

Off	
AlwaysOn	
Gated_Low_Active	
Gated_High_Active	

10.2.1.57 ReferenceElectrodeModeEnumNet enum [ReferenceElectrodeModeEnumNet](#) [strong]

enumerates the electrode subtraction modes

Enumerator

SubtractionOff	
SubtractFromAllOther	
SubtractFromReferenceElectrodeOnly	
SubtractFromAll	

10.2.1.58 ReferenceElectrodeSwitchPositionEnumNet enum [ReferenceElectrodeSwitchPositionEnumNet](#) [strong]

enumerates the possible positions of the reference electrode switch of the ME2100 device

Enumerator

off	
Ref8	
Ref16	
Ref24	
Ref32	

10.2.1.59 RetriggerActionEnumNet enum [RetriggerActionEnumNet](#) [strong]

Enumerates possible retrigger actions for STG200x devices.

Enumerator

raStop	
raRestart	
raIgnore	
raGate	
raSingle	

10.2.1.60 RoboCurrentModeEnumNet enum [RoboCurrentModeEnumNet](#) [strong]

Enumerator

Off	
Break	
Standby	
Reference	
Movement	

10.2.1.61 SampleDstSizeNet enum [SampleDstSizeNet](#) [strong]

Enumerates the destination data format for ChannelBlock functions.

Enumerator

SampleDstSize16	
SampleDstSize32	

10.2.1.62 SampleSizeNet enum [SampleSizeNet](#) [strong]

Enumerates the data format for ChannelBlock functions.

Enumerator

SampleSize16Unsigned	
SampleSize16Signed	
SampleSize24Unsigned	
SampleSize24Signed	
SampleSize32Unsigned	
SampleSize32Signed	
SampleSize64Unsigned	
SampleSize64Signed	

10.2.1.63 SCU_HeadstageIdEnumNet enum [SCU_HeadstageIdEnumNet](#) [strong]

Enumerator

DeviceNotConnected	
Me2100_32Xilinx	
Me2100_32PICiCE40	
Mea2100Mini60PICiCE40	
Mea2100BetaScreen	
Me2100UPA32	
MultiwellMini	
Mea2100Mini120	
Mea2100Mini60ECP5	
eCube	
Me2100Graphene16_32	
GrapheneASIC	
WholeCellPatch	
DeviceHasNoHeadstage	

10.2.1.64 SCUDacqGroupChannelEnumNet enum [SCUDacqGroupChannelEnumNet](#) [strong]

Enumerates the SCU Channel Groups of Datastream

Enumerator

SCU1ElectrodeGroupHS1	
SCU1ElectrodeGroupHS2	
SCU1ElectrodeGroupHS3	
SCU1ElectrodeGroupHS4	
SCU2ElectrodeGroupHS1	
SCU2ElectrodeGroupHS2	
SCU2ElectrodeGroupHS3	
SCU2ElectrodeGroupHS4	
InterfaceADCGroup	
STG1DACSignalGroup	
STG2DACSignalGroup	
DSPAnalogGroup	
DSPDigitalGroup	
IFDigChannelsGroup	
STG1TriggerStatusGroup	
STG1SidebandsGroup	
STG2TriggerStatusGroup	
STG2SidebandsGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.65 SCUDigitalSourceEnumNet enum [SCUDigitalSourceEnumNet](#) [strong]

Enumerates the digital source of the SCU device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
HS1Trigger1Status	

Enumerator

HS1Trigger2Status	
HS1Trigger3Status	
HS1Trigger4Status	
HS1Trigger5Status	
HS1Trigger6Status	
HS1Trigger7Status	
HS1Trigger8Status	
HS1Trigger9Status	
HS1Trigger10Status	
HS1Trigger11Status	
HS1Trigger12Status	
HS1Sideband1	
HS1Sideband2	
HS1Sideband3	
HS1Sideband4	
HS1Sideband5	
HS1Sideband6	
HS1Sideband7	
HS1Sideband8	
HS1Sideband9	
HS1Sideband10	
HS1Sideband11	
HS1Sideband12	
HS2Trigger1Status	
HS2Trigger2Status	
HS2Trigger3Status	
HS2Trigger4Status	
HS2Trigger5Status	
HS2Trigger6Status	
HS2Trigger7Status	
HS2Trigger8Status	
HS2Trigger9Status	
HS2Trigger10Status	
HS2Trigger11Status	
HS2Trigger12Status	
HS2Sideband1	
HS2Sideband2	
HS2Sideband3	
HS2Sideband4	
HS2Sideband5	
HS2Sideband6	
HS2Sideband7	
HS2Sideband8	
HS2Sideband9	
HS2Sideband10	
HS2Sideband11	
HS2Sideband12	
LastPosition	

10.2.1.66 Stg200xDigoutModeEnumNet enum [Stg200xDigoutModeEnumNet](#) [strong]

Enumerates the DigoutMode on STG400x devices.

Enumerator

Monitor	Monitor digital input pins. Digital out is a mirror of the digital input pins.
Manual	Manually set the value on the digital out pins with SetDigoutValue.
SYNCOUT1	show bit 7 to 15 of syncout channel 1 on the digital outputs
SYNCOUT2	show bit 7 to 15 of syncout channel 2 on the digital outputs
SYNCOUT3	show bit 7 to 15 of syncout channel 3 on the digital outputs
SYNCOUT4	show bit 7 to 15 of syncout channel 4 on the digital outputs
SYNCOUT5	show bit 7 to 15 of syncout channel 5 on the digital outputs
SYNCOUT6	show bit 7 to 15 of syncout channel 6 on the digital outputs
SYNCOUT7	show bit 7 to 15 of syncout channel 7 on the digital outputs
SYNCOUT8	show bit 7 to 15 of syncout channel 8 on the digital outputs

10.2.1.67 Stg200xSegmentFlagsEnumNet enum [Stg200xSegmentFlagsEnumNet](#) [strong]

Enumerates Segmentflag options for STG400x devices.

Enumerator

None	No Flags.
UpdateTrigger	Assign all channels to the trigger which number is the given segment number.
DownloadOnly	Only switch the segment for the next download, keep current segment running.
TriggerOnly	Only switch the segment for the next sweep, keep current download segment.
SyncStart	Delay the start the new segment with SendSegmentStart() until the next sweep has finished.

10.2.1.68 Stg200xTriggerStatusEnumNet enum [Stg200xTriggerStatusEnumNet](#) [strong]

Enumerates the STG download mode trigger status

The STG maintains the status for each of the STG200x_NUM_TRIGGER triggers

Enumerator

Idle	
Running	
Finished	
Armed	

10.2.1.69 STG_DestinationEnumNet enum [STG_DestinationEnumNet](#) [strong]

Enumerates the destination for STG downloads.

Enumerator

channeldata_voltage	
channeldata_current	
syncoutdata	
channeldata_positive_voltage	
channeldata_positive_current	
rawdata	
channeldata_current_own_sync	
channeldata_positive_current_own_sync	
channeldata_current_own_boost_gnd_sync	
channeldata_positive_current_own_boost_gnd_sync	
channeldata_current_always_boost	
channeldata_current_always_boost_own_sync	

10.2.1.70 StimulationLayoutConfigurationEnumNet enum [StimulationLayoutConfigurationEnumNet](#) [strong]

enumerates the layout configuration for the MEA2100-256 device

Enumerator

SingleWell	
SixWell	
NineWell	

10.2.1.71 TBSI_DACQDigitalSourceEnumNet enum [TBSI_DACQDigitalSourceEnumNet](#) [strong]

Enumerates the digital source of the TBSI-DACQ device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	

Enumerator

DeviceRunStatus	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
HS1DigitalData1	
HS2DigitalData1	
LastPosition	

10.2.1.72 TcxDeviceTypeEnumNet enum [TcxDeviceTypeEnumNet](#) [strong]

Enumerates the type of TCX devices.

Enumerator

Unknown	
Regular	
BMI	
Nanion	
Warner	

10.2.1.73 TcxSensorTypeEnumNet enum [TcxSensorTypeEnumNet](#) [strong]

Enumerates the sensor types for TCX devices

Enumerator

Reserved5	
Reserved4	
Reserved3	
Reserved2	
Reserved1	
NTC10K	
PT1000	
PT100	

10.2.1.74 TeerClampModeEnumNet enum [TeerClampModeEnumNet](#) [strong]

Enumerator

ClampModeVoltage	
------------------	--

Enumerator

ClampModeCurrent	
ClampModeOpen	
ClampModeInternalCalibration	

10.2.1.75 TeerWaveformEnumNet `enum TeerWaveformEnumNet [strong]`

Enumerator

Rectangle	
Sine	

10.2.1.76 TriggerSourceEnumNet `enum TriggerSourceEnumNet [strong]`

Enumerates the trigger source of the MEA2100 device.

Enumerator

tsNone	
tsDigitalIn1	
tsDigitalIn2	
tsDigitalIn3	
tsDigitalIn4	
tsDigitalIn5	
tsDigitalIn6	
tsDigitalIn7	
tsDigitalIn8	
tsDigitalIn9	
tsDigitalIn10	
tsDigitalIn11	
tsDigitalIn12	
tsDigitalIn13	
tsDigitalIn14	
tsDigitalIn15	
tsDigitalIn16	
tsDigitalIn17	
tsDigitalIn18	
tsDigitalIn19	
tsDigitalIn20	
tsDigitalIn21	
tsDigitalIn22	
tsDigitalIn23	
tsDigitalIn24	

Enumerator

tsDigitalIn25	
tsDigitalIn26	
tsDigitalIn27	
tsDigitalIn28	
tsDigitalIn29	
tsDigitalIn30	
tsDigitalIn31	
tsDigitalIn32	
tsFeedback1	
tsFeedback2	
tsFeedback3	
tsFeedback4	
tsFeedback5	
tsFeedback6	
tsFeedback7	
tsFeedback8	
tsFeedback9	
tsFeedback10	
tsFeedback11	
tsFeedback12	
tsFeedback13	
tsFeedback14	
tsFeedback15	
tsFeedback16	
tsFeedback17	
tsFeedback18	
tsFeedback19	
tsFeedback20	
tsFeedback21	
tsFeedback22	
tsFeedback23	
tsFeedback24	
tsFeedback25	
tsFeedback26	
tsFeedback27	
tsFeedback28	
tsFeedback29	
tsFeedback30	
tsFeedback31	
tsFeedback32	
tsAuxIn1	
tsAuxIn2	
tsDigitalPuse0	
tsDigitalPuse1	
tsDigitalPuse2	
tsDigitalPuse3	
tsDigitalPuse4	
tsDigitalPuse5	
tsDigitalPuse6	
tsDigitalPuse7	

Enumerator

tsDigitalPuse8	
tsDigitalPuse9	
tsDigitalPuse10	
tsDigitalPuse11	
tsDigitalPuse12	
tsDigitalPuse13	
tsDigitalPuse14	
tsDigitalPuse15	
tsDigitalPuse16	
tsDigitalPuse17	
tsDigitalPuse18	
tsDigitalPuse19	
tsDigitalPuse20	
tsDigitalPuse21	
tsDigitalPuse22	
tsDigitalPuse23	
tsDigitalPuse24	
tsDigitalPuse25	
tsDigitalPuse26	
tsDigitalPuse27	
tsDigitalPuse28	
tsDigitalPuse29	
tsDigitalPuse30	
tsDigitalPuse31	
tsTriggered	
tsSidebandBit8	
tsDACQCy1Dev1Runs	
tsDACQCy1Dev2Runs	
tsDACQCy2Dev1Runs	
tsDACQCy2Dev2Runs	

10.2.1.77 UsbVendorIdEnumNet enum `UsbVendorIdEnumNet` [strong]

Enumerator

Unknown	
None	
Renesas	
ASMedia	
Intel	

10.2.1.78 UssingClampModeEnumNet enum `UssingClampModeEnumNet` [strong]

Enumerator

VoltageClamp	
CurrentClamp	
OpenClamp	
Standby	
ElectrodeOffset	

10.2.1.79 UssingUnitEnumNet enum [UssingUnitEnumNet](#) [strong]

Enumerator

Volt	
Ampere	
State	

10.2.1.80 VendorIdEnumNet enum [VendorIdEnumNet](#) [strong]

Enumerates the group of MCS devices to connect to.

Enumerator

Any	
None	
MCS	
PCI	
Cypress	
ALA_VC3	

10.2.1.81 W2100_Accel_Gyro_Select_EnumNet enum [W2100_Accel_Gyro_Select_EnumNet](#) [strong]

enumerates the accelerometer configuration on the W2100 device

Enumerator

Off	
GyroOnly	
AccelOnly	
Both	

10.2.1.82 W2100DacqGroupChannelEnumNet enum [W2100DacqGroupChannelEnumNet](#) [strong]

Enumerates the W2100 Channel Groups of Datastream

Enumerator

InterfaceADCGroup	
DSPDataGroup	
WirelessHeadStageAnalogRE1HS1	
WirelessHeadStageStatusRE1HS1	
WirelessHeadStageAnalogRE1HS2	
WirelessHeadStageStatusRE1HS2	
WirelessHeadStageAnalogRE1HS3	
WirelessHeadStageStatusRE1HS3	
WirelessHeadStageAnalogRE1HS4	
WirelessHeadStageStatusRE1HS4	
WirelessHeadStageAnalogRE2HS1	
WirelessHeadStageStatusRE2HS1	
WirelessHeadStageAnalogRE2HS2	
WirelessHeadStageStatusRE2HS2	
WirelessHeadStageAnalogRE2HS3	
WirelessHeadStageStatusRE2HS3	
WirelessHeadStageAnalogRE2HS4	
WirelessHeadStageStatusRE2HS4	
WirelessHeadStageGyroDataRE1HS1	
WirelessHeadStageAccDataRE1HS1	
WirelessHeadStageGyroDataRE1HS2	
WirelessHeadStageAccDataRE1HS2	
WirelessHeadStageGyroDataRE1HS3	
WirelessHeadStageAccDataRE1HS3	
WirelessHeadStageGyroDataRE1HS4	
WirelessHeadStageAccDataRE1HS4	
WirelessHeadStageGyroDataRE2HS1	
WirelessHeadStageAccDataRE2HS1	
WirelessHeadStageGyroDataRE2HS2	
WirelessHeadStageAccDataRE2HS2	
WirelessHeadStageGyroDataRE2HS3	
WirelessHeadStageAccDataRE2HS3	
WirelessHeadStageGyroDataRE2HS4	
WirelessHeadStageAccDataRE2HS4	
WirelessHeadStageOptoStimCurrentRE1HS1	
WirelessHeadStageReservedARE1HS1	
WirelessHeadStageOptoStimCurrentRE1HS2	
WirelessHeadStageReservedARE1HS2	
WirelessHeadStageOptoStimCurrentRE1HS3	
WirelessHeadStageReservedARE1HS3	
WirelessHeadStageOptoStimCurrentRE1HS4	
WirelessHeadStageReservedARE1HS4	
WirelessHeadStageOptoStimCurrentRE2HS1	

Enumerator

WirelessHeadStageReservedARE2HS1	
WirelessHeadStageOptoStimCurrentRE2HS2	
WirelessHeadStageReservedARE2HS2	
WirelessHeadStageOptoStimCurrentRE2HS3	
WirelessHeadStageReservedARE2HS3	
WirelessHeadStageOptoStimCurrentRE2HS4	
WirelessHeadStageReservedARE2HS4	
WirelessHeadStageReservedBRE1HS1	
WirelessHeadStageReservedCRE1HS1	
WirelessHeadStageReservedBRE1HS2	
WirelessHeadStageReservedCRE1HS2	
WirelessHeadStageReservedBRE1HS3	
WirelessHeadStageReservedCRE1HS3	
WirelessHeadStageReservedBRE1HS4	
WirelessHeadStageReservedCRE1HS4	
WirelessHeadStageReservedBRE2HS1	
WirelessHeadStageReservedCRE2HS1	
WirelessHeadStageReservedBRE2HS2	
WirelessHeadStageReservedCRE2HS2	
WirelessHeadStageReservedBRE2HS3	
WirelessHeadStageReservedCRE2HS3	
WirelessHeadStageReservedBRE2HS4	
WirelessHeadStageReservedCRE2HS4	
IFDigChannelsGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.83 W2100DigitalSourceEnumNet enum W2100DigitalSourceEnumNet [strong]

Enumerates the digital source of the W2100 device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
PulseGenerator	
DigDataFromReceiver	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
DigStreamFromReceiver	
DigStreamPosition	

10.2.1.84 WvcDisplayModeEnumNet enum [WvcDisplayModeEnumNet](#) [strong]

enumerates Wvc display mode

Enumerator

Work	
PC	
Settings	
TouchTest	

10.2.1.85 WvcValveModeEnumNet enum [WvcValveModeEnumNet](#) [strong]

enumerates Wvc valve mode

Enumerator

Manual	
Digital	
Analog	
Table	

10.2.2 Function Documentation

10.2.2.1 OnChannelData() public delegate void Mcs::Usb::OnChannelData (
 [CMcsUsbDacqNet](#)[^] *dacq*,
 int *CbHandle*,
 int *numFrames*)

10.2.2.2 OnDeviceArrivalRemoval() public delegate void Mcs::Usb::OnDeviceArrivalRemoval (
 [CMcsUsbListEntryNet](#)[^] *entry*)

Delegate to show a device arrival or removal.

10.2.2.3 OnError() public delegate void Mcs::Usb::OnError (
 String[^] *msg*,
 int *action*)

10.2.2.4 OnMcsUsbDeviceState() public delegate void OnMcsUsbDeviceState (
usbSetupPacket_t^ request)

10.2.2.5 OnMcsUsbDeviceStateCallback() private delegate void OnMcsUsbDeviceStateCallback (
IntPtr pThis,
uint32_t size,
IntPtr buffer)

10.2.2.6 OnMwPollStatus() public delegate void Mcs::Usb::OnMwPollStatus (
unsigned int CurrentTemp,
unsigned int PlateState,
unsigned int SwitchState)

10.2.2.7 OnStg200xDataHandler() public delegate void Mcs::Usb::OnStg200xDataHandler (
uint32_t trigger)

10.2.2.8 OnStg200xErrorHandler() public delegate void Mcs::Usb::OnStg200xErrorHandler ()

10.2.2.9 OnStgPollStatus() public delegate void Mcs::Usb::OnStgPollStatus (
unsigned int status,
StgStatusNet^ stgStatusNet,
array< int >^ index_list)

10.2.2.10 OnUpdateFirmwareProgress() public delegate void Mcs::Usb::OnUpdateFirmwareProgress
(
int)

10.2.2.11 OnUpdateFirmwareStatusChange() public delegate void Mcs::Usb::OnUpdateFirmware←
StatusChange (
String^)

10.2.2.12 RoboStatusEventDelegate() public delegate void Mcs::Usb::RoboStatusEventDelegate (
array< unsigned char >^ buffer)

11 Class Documentation

11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference

Public Attributes

- [W2100DacqGroupChannelEnumNet](#) `dacqgroup`
- `int` `channel`
- `int` `amplification`

11.1.1 Member Data Documentation

11.1.1.1 `amplification` `int` `amplification`

11.1.1.2 `channel` `int` `channel`

11.1.1.3 `dacqgroup` [W2100DacqGroupChannelEnumNet](#) `dacqgroup`

11.2 BatteryState Class Reference

Properties

- `unsigned int` [Charge](#) `[get]`
- `unsigned int` [Voltage](#) `[get]`
- `System::String^` [ChargeString](#) `[get]`
- `System::String^` [ChargeRegionString](#) `[get]`
- `System::String^` [VoltageString](#) `[get]`

11.2.1 Property Documentation

11.2.1.1 `Charge` `unsigned int` `Charge` `[get]`

11.2.1.2 `ChargeRegionString` `System::String^` `ChargeRegionString` `[get]`


11.2.1.3 ChargeString `System:: String^ ChargeString [get]`

11.2.1.4 Voltage `unsigned int Voltage [get]`

11.2.1.5 VoltageString `System:: String^ VoltageString [get]`

11.3 BesselFilterHighPassNet Class Reference

Inheritance diagram for BesselFilterHighPassNet:

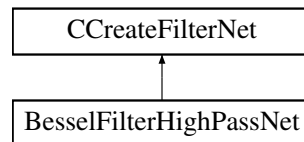

Public Member Functions

- [BesselFilterHighPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)


Additional Inherited Members

11.3.1 Constructor & Destructor Documentation

11.3.1.1 BesselFilterHighPassNet() `BesselFilterHighPassNet (`
 int numCoefSets,
 int order,
 double sampleRate,
 double cutoffFrequency,
 double scale)

11.4 BesselFilterLowPassNet Class Reference

Inheritance diagram for BesselFilterLowPassNet:

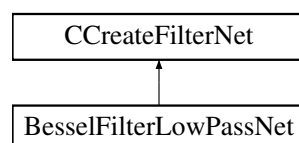

Public Member Functions

- [BesselFilterLowPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)


Additional Inherited Members

11.4.1 Constructor & Destructor Documentation

11.4.1.1 BesselFilterLowPassNet() `BesselFilterLowPassNet (`
 int numCoefSets,
 int order,
 double sampleRate,
 double cutoffFrequency,
 double scale)

11.5 ButterworthFilterHighPassNet Class Reference

Inheritance diagram for ButterworthFilterHighPassNet:

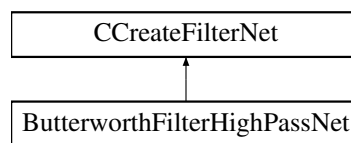

Public Member Functions

- [ButterworthFilterHighPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)


Additional Inherited Members

11.5.1 Constructor & Destructor Documentation

11.5.1.1 ButterworthFilterHighPassNet() `ButterworthFilterHighPassNet (`
 int numCoefSets,
 int order,
 double sampleRate,
 double cutoffFrequency,
 double scale)

11.6 ButterworthFilterLowPassNet Class Reference

Inheritance diagram for ButterworthFilterLowPassNet:

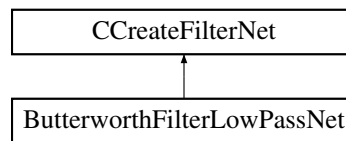

Public Member Functions

- [ButterworthFilterLowPassNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members

11.6.1 Constructor & Destructor Documentation

11.6.1.1 ButterworthFilterLowPassNet() [ButterworthFilterLowPassNet](#) (


```

    int numCoefSets,
    int order,
    double sampleRate,
    double cutoffFrequency,
    double scale )

```

11.7 CChannelTestDeviceNet Class Reference

Inheritance diagram for CChannelTestDeviceNet:

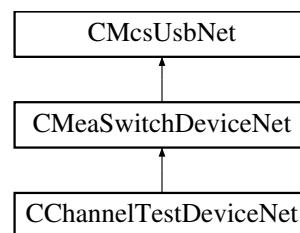

Public Member Functions

- [CChannelTestDeviceNet](#) ()
- [~CChannelTestDeviceNet](#) ()
- void [SetWaveform](#) (unsigned int Waveform)
- void [SetAmplitude](#) (unsigned int Amplitude)
- void [SetFrequency](#) (unsigned int Frequency)
- void [SetAttenuation](#) (unsigned int Attenuation)

11.7.1.1 CChannelTestDeviceNet() CChannelTestDeviceNet ()

11.7.1.2 ~CChannelTestDeviceNet() ~CChannelTestDeviceNet ()

```
11.7.2.1 SetAmplitude() void SetAmplitude (
    unsigned int Amplitude )
```

```
11.7.2.2 SetAttenuation() void SetAttenuation (
    unsigned int Attenuation )
```

11.7.2.3 SetFrequency()

```
void SetFrequency (
    unsigned int Frequency )
```

```
11.7.2.4 SetWaveform() void SetWaveform (
    unsigned int Waveform )
```

11.8 CCMOSMea_FunctionNet Class Reference

```

graph BT
    CCMOSMea_FunctionNet --> CMcsUsbFunctionNet

```

Public Member Functions

- [CCMOSMea_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ cMOSMea_↵
FunctionPointerContainer)
- [CCMOSMea_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- void [SetADCInputOffset](#) (int32_t offset)
- int32_t [GetADCInputOffset](#) ()
- void [SetSourceDrain](#) (int32_t voltage)
- int32_t [GetSourceDrain](#) ()
- void [SetSourceGate](#) (int32_t voltage)
- int32_t [GetSourceGate](#) ()
- void [SetSourceBulk](#) (int32_t voltage)
- int32_t [GetSourceBulk](#) ()
- void [SetGate](#) (int32_t voltage)
- int32_t [GetGate](#) ()
- void [SetBath](#) (int32_t voltage)
- int32_t [GetBath](#) ()
- int32_t [GetGNDI](#) ()
- int32_t [GetVDDI](#) ()
- int32_t [GetVDD3I](#) ()
- void [UpdateTransistorVoltages](#) ()
- bool [AreTransistorVoltagesSet](#) ()
- void [PowerChip](#) (bool on)
- bool [IsChipPowered](#) ()
- [enCMosMeaChipType](#) [DetectChipType](#) ()
- void [SetGateToVOP](#) ()
- void [SetGateFloating](#) ()
- bool [IsGateFloating](#) ()
- void [VOPSTimerSetResetTimes](#) (uint32_t ResetTime, uint32_t IntervalTime)
- void [VOPSTimerSetResetTimes](#) (uint32_t ResetTime, uint32_t IntervalTime, uint32_t HPFilterResetTime)
- void [SetBathMode](#) ([CMOSMeaBathModeEnumNet](#) Mode)
- [CMOSMeaBathModeEnumNet](#) [GetBathMode](#) ()
- void [SetNeurochipMemoryData](#) (uint16_t MemAddress, uint32_t MemData)
- void [SetNeurochipMemoryData](#) (uint16_t MemAddress, array< uint32_t >^ MemData)
- uint32_t [GetNeurochipMemoryData](#) (uint16_t MemAddress)
- array< uint32_t > ^ [GetNeurochipMemoryData](#) (uint16_t MemAddress, uint32_t RequestLength)
- uint32_t [GetNeurochipMemorySize](#) ()
- uint32_t [GetMaxNumOfColumns](#) (uint32_t Samplerate)
- void [SetStimulusSites](#) (List< int16_t >^ SwitchPosition)
- List< int16_t > ^ [GetStimulusSites](#) ()
- void [ClearSTGOutput](#) (uint32_t Channel)
- uint32_t [GetNumberOfSupportedGroups](#) ()
- uint32_t [GetNumberOfSupportedGroups](#) (uint32_t virtualDevice)
- [DacqGroupChannelEnumNet](#) [GetGroupID](#) (uint32_t Index)
- [DacqGroupChannelEnumNet](#) [GetGroupID](#) (uint32_t Index, uint32_t virtualDevice)
- uint32_t [GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumNet](#) GroupID)
- uint32_t [GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumNet](#) GroupID, uint32_t virtualDevice)
- [DacqMeaGroupTypeEnumNet](#) [GetGroupType](#) ([DacqGroupChannelEnumNet](#) GroupID)
- [DacqMeaGroupTypeEnumNet](#) [GetGroupType](#) ([DacqGroupChannelEnumNet](#) GroupID, uint32_t virtual↵
Device)
- void [EnableChannelsInGroup](#) ([DacqGroupChannelEnumNet](#) GroupID, List< bool >^ EnabledChannelsBit↵
Map)
- void [EnableChannelsInGroup](#) ([DacqGroupChannelEnumNet](#) GroupID, List< bool >^ EnabledChannelsBit↵
Map, uint32_t virtualDevice)
- List< bool > ^ [GetEnabledChannelsInGroup](#) ([DacqGroupChannelEnumNet](#) GroupID)

- `List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID)`
- `SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `uint32_t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID)`
- `uint32_t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID)`
- `CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `int32_t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID)`
- `int32_t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `int32_t GetGroupADCBits (DacqGroupChannelEnumNet GroupID)`
- `int32_t GetGroupADCBits (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)`
- `uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel↵ Number)`
- `uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel↵ Number, uint32_t virtualDevice)`
- `CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber)`
- `CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaIFDigChannelEnumNet GetGroupChannelBitmaskIFDigChannels (uint32_t ChannelNumber)`
- `CMOSMeaIFDigChannelEnumNet GetGroupChannelBitmaskIFDigChannels (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber)`
- `CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber)`
- `CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel↵ Number)`
- `CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel↵ Number, uint32_t virtualDevice)`
- `CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t ChannelNumber)`
- `CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32_t ChannelNumber)`
- `CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32_t ChannelNumber, uint32_t virtualDevice)`
- `CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel↵ Number)`
- `CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel↵ Number, uint32_t virtualDevice)`
- `CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_↵ t ChannelNumber)`
- `CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_↵ t ChannelNumber, uint32_t virtualDevice)`

Additional Inherited Members

11.8.1 Constructor & Destructor Documentation

11.8.1.1 CCMOSMea_FunctionNet() [1/2] `CCMOSMea_FunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ cMOSMea_FunctionPointerContainer)`

11.8.1.2 CCMOSMea_FunctionNet() [2/2] `CCMOSMea_FunctionNet (CMcsUsbNet^ mcsusb)`

11.8.2 Member Function Documentation

11.8.2.1 AreTransistorVoltagesSet() `bool AreTransistorVoltagesSet ()`

11.8.2.2 ClearSTGOutput() `void ClearSTGOutput (uint32_t Channel)`

11.8.2.3 DetectChipType() `enCMosMeaChipType DetectChipType ()`

11.8.2.4 EnableChannelsInGroup() [1/2] `void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBitMap)`

11.8.2.5 EnableChannelsInGroup() [2/2] `void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBitMap, uint32_t virtualDevice)`

11.8.2.6 GetADCInputOffset() `int32_t GetADCInputOffset ()`

11.8.2.7 GetBath() `int32_t GetBath ()`

11.8.2.8 GetBathMode() `CMOSMeaBathModeEnumNet` GetBathMode ()

11.8.2.9 GetEnabledChannelsInGroup() [1/2] `List<bool>` ^ GetEnabledChannelsInGroup (
`DacqGroupChannelEnumNet` *GroupID*)

11.8.2.10 GetEnabledChannelsInGroup() [2/2] `List<bool>` ^ GetEnabledChannelsInGroup (
`DacqGroupChannelEnumNet` *GroupID*,
`uint32_t` *virtualDevice*)

11.8.2.11 GetGate() `int32_t` GetGate ()

11.8.2.12 GetGNDI() `int32_t` GetGNDI ()

11.8.2.13 GetGroupADCBits() [1/2] `int32_t` GetGroupADCBits (
`DacqGroupChannelEnumNet` *GroupID*)

11.8.2.14 GetGroupADCBits() [2/2] `int32_t` GetGroupADCBits (
`DacqGroupChannelEnumNet` *GroupID*,
`uint32_t` *virtualDevice*)

11.8.2.15 GetGroupChannelBitmaskBySelect() [1/2] `uint32_t` GetGroupChannelBitmaskBySelect (
`DacqGroupChannelEnumNet` *GroupID*,
`uint32_t` *ChannelNumber*)

11.8.2.16 GetGroupChannelBitmaskBySelect() [2/2] `uint32_t` GetGroupChannelBitmaskBySelect (
`DacqGroupChannelEnumNet` *GroupID*,
`uint32_t` *ChannelNumber*,
`uint32_t` *virtualDevice*)

11.8.2.17 GetGroupChannelBitmaskHS1NCBathCurrent() [1/2] [CMOSMeaHeadstage1NCBathCurrentEnumNet](#)

```
GetGroupChannelBitmaskHS1NCBathCurrent (
    uint32_t ChannelNumber )
```

11.8.2.18 GetGroupChannelBitmaskHS1NCBathCurrent() [2/2] [CMOSMeaHeadstage1NCBathCurrentEnumNet](#)

```
GetGroupChannelBitmaskHS1NCBathCurrent (
    uint32_t ChannelNumber,
    uint32_t virtualDevice )
```

11.8.2.19 GetGroupChannelBitmaskHS1NCCol2Current() [1/2] [CMOSMeaHeadstage1NCCol2CurrentEnumNet](#)

```
GetGroupChannelBitmaskHS1NCCol2Current (
    uint32_t ChannelNumber )
```

11.8.2.20 GetGroupChannelBitmaskHS1NCCol2Current() [2/2] [CMOSMeaHeadstage1NCCol2CurrentEnumNet](#)

```
GetGroupChannelBitmaskHS1NCCol2Current (
    uint32_t ChannelNumber,
    uint32_t virtualDevice )
```

11.8.2.21 GetGroupChannelBitmaskHS1NChipTemp() [1/2] [CMOSMeaHeadstage1NChipTempEnumNet](#) Get↔

```
GroupChannelBitmaskHS1NChipTemp (
    uint32_t ChannelNumber )
```

11.8.2.22 GetGroupChannelBitmaskHS1NChipTemp() [2/2] [CMOSMeaHeadstage1NChipTempEnumNet](#) Get↔

```
GroupChannelBitmaskHS1NChipTemp (
    uint32_t ChannelNumber,
    uint32_t virtualDevice )
```

11.8.2.23 GetGroupChannelBitmaskHS1Sidebands() [1/2] [CMOSMeaHS1SidebandEnumNet](#) GetGroup↔

```
ChannelBitmaskHS1Sidebands (
    uint32_t ChannelNumber )
```

11.8.2.24 GetGroupChannelBitmaskHS1Sidebands() [2/2] [CMOSMeaHS1SidebandEnumNet](#) GetGroup↔

```
ChannelBitmaskHS1Sidebands (
    uint32_t ChannelNumber,
    uint32_t virtualDevice )
```

11.8.2.25 GetGroupChannelBitmaskHS1TriggerStatus() [1/2] [CMOSMeaHS1TriggerStatusEnumNet](#) Get↔
GroupChannelBitmaskHS1TriggerStatus (
 uint32_t ChannelNumber)

11.8.2.26 GetGroupChannelBitmaskHS1TriggerStatus() [2/2] [CMOSMeaHS1TriggerStatusEnumNet](#) Get↔
GroupChannelBitmaskHS1TriggerStatus (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

11.8.2.27 GetGroupChannelBitmaskIFDigChannels() [1/2] [CMOSMeaIFDigChannelEnumNet](#) GetGroup↔
ChannelBitmaskIFDigChannels (
 uint32_t ChannelNumber)

11.8.2.28 GetGroupChannelBitmaskIFDigChannels() [2/2] [CMOSMeaIFDigChannelEnumNet](#) GetGroup↔
ChannelBitmaskIFDigChannels (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

11.8.2.29 GetGroupChannelBitmaskInterfaceADC() [1/2] [CMOSMeaInterfaceADCEnumNet](#) GetGroup↔
ChannelBitmaskInterfaceADC (
 uint32_t ChannelNumber)

11.8.2.30 GetGroupChannelBitmaskInterfaceADC() [2/2] [CMOSMeaInterfaceADCEnumNet](#) GetGroup↔
ChannelBitmaskInterfaceADC (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

11.8.2.31 GetGroupChannelBitmaskPacketFrameContext() [1/2] [CMOSMeaPacketFrameContextGroupEnumNet](#)
GetGroupChannelBitmaskPacketFrameContext (
 uint32_t ChannelNumber)

11.8.2.32 GetGroupChannelBitmaskPacketFrameContext() [2/2] [CMOSMeaPacketFrameContextGroupEnumNet](#)
GetGroupChannelBitmaskPacketFrameContext (
 uint32_t ChannelNumber,
 uint32_t virtualDevice)

11.8.2.33 GetGroupChannelBitmaskSTG1DACSignal() [1/2] [CMOSMeaSTG1DACSignalEnumNet](#) GetGroup↔
 ChannelBitmaskSTG1DACSignal (
 uint32_t *ChannelNumber*)

11.8.2.34 GetGroupChannelBitmaskSTG1DACSignal() [2/2] [CMOSMeaSTG1DACSignalEnumNet](#) GetGroup↔
 ChannelBitmaskSTG1DACSignal (
 uint32_t *ChannelNumber*,
 uint32_t *virtualDevice*)

11.8.2.35 GetGroupDCOffset() [1/2] int32_t GetGroupDCOffset (
 [DacqGroupChannelEnumNet](#) *GroupID*)

11.8.2.36 GetGroupDCOffset() [2/2] int32_t GetGroupDCOffset (
 [DacqGroupChannelEnumNet](#) *GroupID*,
 uint32_t *virtualDevice*)

11.8.2.37 GetGroupID() [1/2] [DacqGroupChannelEnumNet](#) GetGroupID (
 uint32_t *Index*)

11.8.2.38 GetGroupID() [2/2] [DacqGroupChannelEnumNet](#) GetGroupID (
 uint32_t *Index*,
 uint32_t *virtualDevice*)

11.8.2.39 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
 [DacqGroupChannelEnumNet](#) *GroupID*)

11.8.2.40 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
 [DacqGroupChannelEnumNet](#) *GroupID*,
 uint32_t *virtualDevice*)

11.8.2.41 GetGroupResolutionPerDigit() [1/2] uint32_t GetGroupResolutionPerDigit (
 [DacqGroupChannelEnumNet](#) *GroupID*)

11.8.2.42 GetGroupResolutionPerDigit() [2/2] `uint32_t GetGroupResolutionPerDigit (`
 `DacqGroupChannelEnumNet GroupID,`
 `uint32_t virtualDevice)`

11.8.2.43 GetGroupSampleSize() [1/2] `SampleSizeNet GetGroupSampleSize (`
 `DacqGroupChannelEnumNet GroupID)`

11.8.2.44 GetGroupSampleSize() [2/2] `SampleSizeNet GetGroupSampleSize (`
 `DacqGroupChannelEnumNet GroupID,`
 `uint32_t virtualDevice)`

11.8.2.45 GetGroupType() [1/2] `DacqMeaGroupTypeEnumNet GetGroupType (`
 `DacqGroupChannelEnumNet GroupID)`

11.8.2.46 GetGroupType() [2/2] `DacqMeaGroupTypeEnumNet GetGroupType (`
 `DacqGroupChannelEnumNet GroupID,`
 `uint32_t virtualDevice)`

11.8.2.47 GetGroupUnit() [1/2] `CMOSMeaValueUnitEnumNet GetGroupUnit (`
 `DacqGroupChannelEnumNet GroupID)`

11.8.2.48 GetGroupUnit() [2/2] `CMOSMeaValueUnitEnumNet GetGroupUnit (`
 `DacqGroupChannelEnumNet GroupID,`
 `uint32_t virtualDevice)`

11.8.2.49 GetMaxNumOfColumns() `uint32_t GetMaxNumOfColumns (`
 `uint32_t Samplerate)`

11.8.2.50 GetNeurochipMemoryData() [1/2] `uint32_t GetNeurochipMemoryData (`
 `uint16_t MemAddress)`

11.8.2.51 GetNeurochipMemoryData() [2/2] `array<uint32_t> ^ GetNeurochipMemoryData (`
 `uint16_t MemAddress,`
 `uint32_t RequestLength)`

11.8.2.52 GetNeurochipMemorySize() `uint32_t GetNeurochipMemorySize ()`

11.8.2.53 GetNumberOfSupportedGroups() [1/2] `uint32_t GetNumberOfSupportedGroups ()`

11.8.2.54 GetNumberOfSupportedGroups() [2/2] `uint32_t GetNumberOfSupportedGroups (`
 `uint32_t virtualDevice)`

11.8.2.55 GetSourceBulk() `int32_t GetSourceBulk ()`

11.8.2.56 GetSourceDrain() `int32_t GetSourceDrain ()`

11.8.2.57 GetSourceGate() `int32_t GetSourceGate ()`

11.8.2.58 GetStimulusSites() `List<int16_t> ^ GetStimulusSites ()`

11.8.2.59 GetVDD3I() `int32_t GetVDD3I ()`

11.8.2.60 GetVDDI() `int32_t GetVDDI ()`

11.8.2.61 IsChipPowered() `bool IsChipPowered ()`

11.8.2.62 IsGateFloating() `bool IsGateFloating ()`

11.8.2.63 PowerChip() `void PowerChip (`
`bool on)`

11.8.2.64 SetADCInputOffset() `void SetADCInputOffset (`
`int32_t offset)`

11.8.2.65 SetBath() `void SetBath (`
`int32_t voltage)`

11.8.2.66 SetBathMode() `void SetBathMode (`
`CMOSMeaBathModeEnumNet Mode)`

11.8.2.67 SetGate() `void SetGate (`
`int32_t voltage)`

11.8.2.68 SetGateFloating() `void SetGateFloating ()`

11.8.2.69 SetGateToVOP() `void SetGateToVOP ()`

11.8.2.70 SetNeurochipMemoryData() [1/2] `void SetNeurochipMemoryData (`
`uint16_t MemAddress,`
`array< uint32_t >^ MemData)`

11.8.2.71 SetNeurochipMemoryData() [2/2] `void SetNeurochipMemoryData (`
`uint16_t MemAddress,`
`uint32_t MemData)`

11.8.2.72 SetSourceBulk() `void SetSourceBulk (`
`int32_t voltage)`

11.8.2.73 SetSourceDrain() `void SetSourceDrain (`
`int32_t voltage)`

11.8.2.74 SetSourceGate() `void SetSourceGate (`
`int32_t voltage)`

11.8.2.75 SetStimulusSites() `void SetStimulusSites (`
`List< int16_t >^ SwitchPosition)`


11.8.2.76 UpdateTransistorVoltages() `void UpdateTransistorVoltages ()`

11.8.2.77 VOPSTimerSetResetTimes() [1/2] `void VOPSTimerSetResetTimes (`
`uint32_t ResetTime,`
`uint32_t IntervalTime)`

11.8.2.78 VOPSTimerSetResetTimes() [2/2] `void VOPSTimerSetResetTimes (`
`uint32_t ResetTime,`
`uint32_t IntervalTime,`
`uint32_t HPFilterResetTime)`

11.9 CCMOSMeaDeviceNet Class Reference

Inheritance diagram for CCMOSMeaDeviceNet:

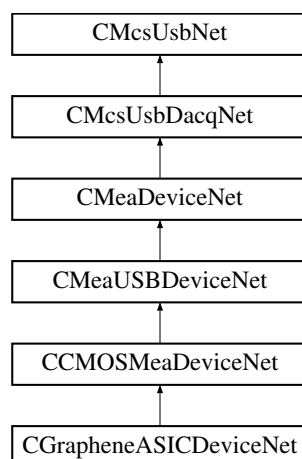

Classes

- class [CRegionOfInterestRect](#)

Public Member Functions

- [CCMOSMeaDeviceNet](#) (void)
- [~CCMOSMeaDeviceNet](#) ()
- virtual void [SetBaseSamplerate](#) (int BaseSamplerate)
- int [GetBaseSamplerate](#) ()
- virtual array< int > ^ [GetAvailableBaseSamplerates](#) ()
- int [GetMaxReadableColumns](#) ()
- virtual void [SetRegionOfInterests](#) (System::Collections::Generic::Dictionary< int, [CRegionOfInterestRect](#) ^>^ rois)
- void [UpdateChannelBlock](#) (int queuesize, int threshold, int channels_in_block)
- System::Collections::Generic::Dictionary< int, array< array< int16_t > ^> ^> ^ [GetCMOSDataDictionary](#) (int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint16_t > ^> ^ [GetChannelDataUI16](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int16_t > ^> ^ [GetChannelDataI16](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint32_t > ^> ^ [GetChannelDataUI32](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int32_t > ^> ^ [GetChannelDataI32](#) (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Properties

- [CCMOSMea_FunctionNet](#) ^ [CMosMea](#) [get]
- [CStimulusFunctionNet](#) ^ [Stimulus](#) [get]

Additional Inherited Members

11.9.1 Constructor & Destructor Documentation

11.9.1.1 CCMOSMeaDeviceNet() [CCMOSMeaDeviceNet](#) (void)

11.9.1.2 ~CCMOSMeaDeviceNet() [~CCMOSMeaDeviceNet](#) ()

11.9.2 Member Function Documentation

11.9.2.1 GetAvailableBaseSamplerates() virtual array<int> ^ GetAvailableBaseSamplerates ()
[virtual]

Reimplemented in [CGrapheneASICDeviceNet](#).

11.9.2.2 GetBaseSamplerate() int GetBaseSamplerate ()

11.9.2.3 GetChannelDataI16() System::Collections::Generic::Dictionary<int, array<int16_t>^> ^
GetChannelDataI16 (
 [DacqGroupChannelEnumNet](#) group,
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret)

11.9.2.4 GetChannelDataI32() System::Collections::Generic::Dictionary<int, array<int32_t>^> ^
GetChannelDataI32 (
 [DacqGroupChannelEnumNet](#) group,
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret)

11.9.2.5 GetChannelDataUI16() System::Collections::Generic::Dictionary<int, array<uint16_t>^> ^
^ GetChannelDataUI16 (
 [DacqGroupChannelEnumNet](#) group,
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret)

11.9.2.6 GetChannelDataUI32() System::Collections::Generic::Dictionary<int, array<uint32_t>^> ^
^ GetChannelDataUI32 (
 [DacqGroupChannelEnumNet](#) group,
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret)

11.9.2.7 GetCMOSDataDictionary() System::Collections::Generic::Dictionary<int, array<array<int16_t>^>^> ^
GetCMOSDataDictionary (
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret)

11.9.2.8 GetMaxReadableColumns() `int GetMaxReadableColumns ()`

11.9.2.9 SetBaseSamplerate() `virtual void SetBaseSamplerate (
int BaseSamplerate) [virtual]`

Reimplemented in [CGrapheneASICDeviceNet](#).

11.9.2.10 SetRegionOfInterests() `virtual void SetRegionOfInterests (
System::Collections::Generic::Dictionary< int, CRegionOfInterestRect>^ rois)
[virtual]`

11.9.2.11 UpdateChannelBlock() `void UpdateChannelBlock (
int queuesize,
int threshold,
int channels_in_block)`


11.9.3 Property Documentation

11.9.3.1 CMosMea [CCMOSMea_FunctionNet](#)^ CMosMea [get]

11.9.3.2 Stimulus [CStimulusFunctionNet](#)^ Stimulus [get]

11.10 CCreateFilterNet Class Reference

Inheritance diagram for CCreateFilterNet:

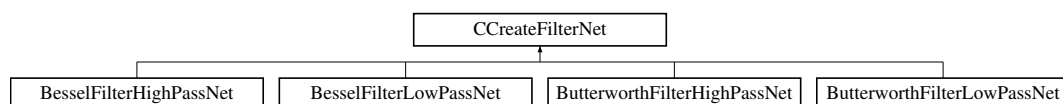

Public Member Functions

- [CCreateFilterNet](#) (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)
- [~CCreateFilterNet](#) ()
- [CFilterCoefficientsNet](#) ^ [GetBiQuad](#) (int index)
- array< [CFilterCoefficientsNet](#) ^> ^ [GetBiQuads](#) ()

Static Public Member Functions

- static int [FindFilter](#) (array< [CFilterCoefficientsNet](#)>^ coef, array< [CCreateFilterNet](#)>^ param)
- static int [FindFilter](#) (array< array< uint64_t >^>^ coef, array< [CCreateFilterNet](#)>^ param, [CFilterCoefficientsNet::s_FilterAttributesNet](#)& FiltAttr, bool DoMCSLegacyCompare)

Protected Member Functions

- [CCreateFilterNet](#) (int numCoefSets, CCreateFilter *pCreateFilter)

Properties

- int [NumCoefSets](#) [get]
- int [Order](#) [get]
- double [SampleRate](#) [get]
- double [CutoffFrequency](#) [get]
- double [Scale](#) [get]

11.10.1 Constructor & Destructor Documentation

11.10.1.1 CCreateFilterNet() [1/2] [CCreateFilterNet](#) (
int *numCoefSets*,
int *order*,
double *sampleRate*,
double *cutoffFrequency*,
double *scale*)

11.10.1.2 ~CCreateFilterNet() [~CCreateFilterNet](#) ()

11.10.1.3 CCreateFilterNet() [2/2] [CCreateFilterNet](#) (
int *numCoefSets*,
CCreateFilter * *pCreateFilter*) [protected]

11.10.2 Member Function Documentation

11.10.2.1 FindFilter() [1/2] static int [FindFilter](#) (
array< array< uint64_t >^>^ *coef*,
array< [CCreateFilterNet](#)>^ *param*,
[CFilterCoefficientsNet::s_FilterAttributesNet](#)& *FiltAttr*,
bool *DoMCSLegacyCompare*) [static]

11.10.2.2 FindFilter() [2/2] static int FindFilter (
array< CFilterCoefficientsNet^>^ coef,
array< CCreateFilterNet^>^ param) [static]

11.10.2.3 GetBiQuad() CFilterCoefficientsNet ^ GetBiQuad (
int index)

11.10.2.4 GetBiQuads() array<CFilterCoefficientsNet^> ^ GetBiQuads ()

11.10.3 Property Documentation

11.10.3.1 CutoffFrequency double CutoffFrequency [get]

11.10.3.2 NumCoefSets int NumCoefSets [get]


11.10.3.3 Order int Order [get]

11.10.3.4 SampleRate double SampleRate [get]

11.10.3.5 Scale double Scale [get]

11.11 CDacCalibrationFunctionNet Class Reference

Inheritance diagram for CDacCalibrationFunctionNet:

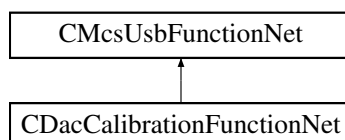

Public Member Functions

- [CDacCalibrationFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pDacCalibrationFunctionPointerContainer)
Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.
- [CDacCalibrationFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CDacCalibrationFunctionNet](#) ()
- [!CDacCalibrationFunctionNet](#) ()
- void [SetDacOffset](#) (uint16_t dacChannel, int32_t offset)
Sets the offset of a DAC channel.
- int32_t [GetDacOffset](#) (uint16_t dacChannel)
Gets the offset of a DAC channel.
- void [BurnDacOffset](#) (uint16_t dacChannel)
Writes the offset of a DAC channel to permanent memory.

Additional Inherited Members**11.11.1 Detailed Description****11.11.2 Constructor & Destructor Documentation**

11.11.2.1 [CDacCalibrationFunctionNet](#)() [1/2] [CDacCalibrationFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pDacCalibrationFunctionPointerContainer)

Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.

11.11.2.2 [CDacCalibrationFunctionNet](#)() [2/2] [CDacCalibrationFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

11.11.2.3 [~CDacCalibrationFunctionNet](#)() virtual [~CDacCalibrationFunctionNet](#) () [virtual]

11.11.2.4 [!CDacCalibrationFunctionNet](#)() [!CDacCalibrationFunctionNet](#) ()

11.11.3 Member Function Documentation

11.11.3.1 [BurnDacOffset](#)() void [BurnDacOffset](#) (
uint16_t dacChannel)

Writes the offset of a DAC channel to permanent memory.

Parameters

<i>dacChannel</i>	The DAC channel number.
-------------------	-------------------------

11.11.3.2 GetDacOffset() `int32_t GetDacOffset (`
 `uint16_t dacChannel)`

Gets the offset of a DAC channel.

Parameters

<i>dacChannel</i>	The DAC channel number.
-------------------	-------------------------

Returns

The offset in digits.

11.11.3.3 SetDacOffset() `void SetDacOffset (`
 `uint16_t dacChannel,`
 `int32_t offset)`


Sets the offset of a DAC channel.

Parameters

<i>dacChannel</i>	The DAC channel number.
<i>offset</i>	The offset in digits.

11.12 CDacqGroupChannelGenericSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelGenericSelectionNet:

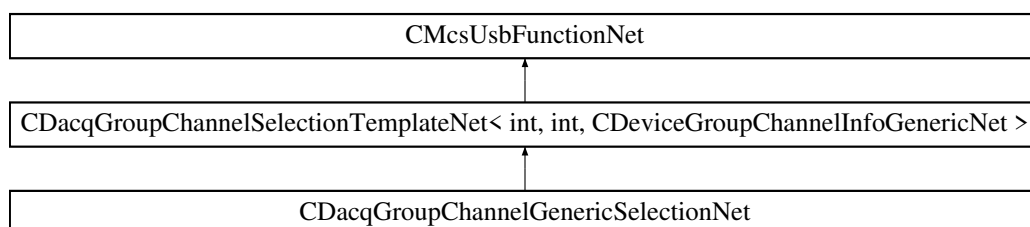**Public Member Functions**

- [CDacqGroupChannelGenericSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)


Additional Inherited Members

11.12.1 Constructor & Destructor Documentation

11.12.1.1 CDacqGroupChannelGenericSelectionNet() `CDacqGroupChannelGenericSelectionNet` (`CMcsUsbNet^ mcsusb`)

11.13 CDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelSelectionNet:

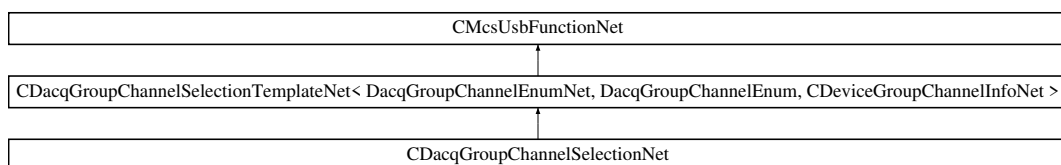

Public Member Functions

- `CDacqGroupChannelSelectionNet` (`CMcsUsbNet^ mcsusb`)


Additional Inherited Members

11.13.1 Constructor & Destructor Documentation

11.13.1.1 CDacqGroupChannelSelectionNet() `CDacqGroupChannelSelectionNet` (`CMcsUsbNet^ mcsusb`)

11.14 CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference

Inheritance diagram for CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >:

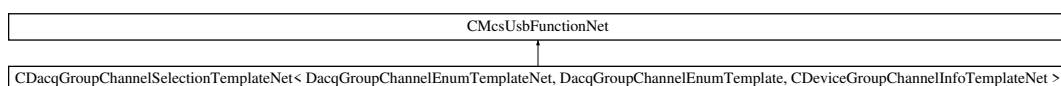

Public Member Functions

- [CDacqGroupChannelSelectionTemplateNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- [uint32_t GetNumberOfSupportedGroups](#) ()
- [uint32_t GetNumberOfSupportedGroups](#) ([uint32_t](#) virtualDevice)
- [DacqGroupChannelEnumTemplateNet GetGroupID](#) ([uint32_t](#) Index)
- [DacqGroupChannelEnumTemplateNet GetGroupID](#) ([uint32_t](#) Index, [uint32_t](#) virtualDevice)
- [uint32_t GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [uint32_t GetGroupNumberOfChannels](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [DacqMeaGroupTypeEnumNet GetGroupType](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [DacqMeaGroupTypeEnumNet GetGroupType](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap)
- [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap, [uint32_t](#) virtualDevice)
- [List< bool >](#)[^] [GetEnabledChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [List< bool >](#)[^] [GetEnabledChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [SampleSizeNet GetGroupSampleSize](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID)
- [SampleSizeNet GetGroupSampleSize](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [uint32_t](#) virtualDevice)
- [List< CDeviceGroupChannelInfoTemplateNet](#)[^][>] [GetDeviceGroupChannelInfos](#) ()
- [List< CDeviceGroupChannelInfoTemplateNet](#)[^][>] [GetDeviceGroupChannelInfos](#) ([uint32_t](#) virtualDevice)

Additional Inherited Members

11.14.1 Constructor & Destructor Documentation

11.14.1.1 CDacqGroupChannelSelectionTemplateNet() [CDacqGroupChannelSelectionTemplateNet](#) ([CMcsUsbNet](#)[^] mcsusb)

11.14.2 Member Function Documentation

11.14.2.1 EnableChannelsInGroup() [1/2] [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap)

11.14.2.2 EnableChannelsInGroup() [2/2] [void EnableChannelsInGroup](#) ([DacqGroupChannelEnumTemplateNet](#) GroupID, [List< bool >](#)[^] EnabledChannelsBitMap, [uint32_t](#) virtualDevice)

11.14 CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference#17

11.14.2.3 GetDeviceGroupChannelInfos() [1/2] List<CDeviceGroupChannelInfoTemplateNet^> ^ Get↔
DeviceGroupChannelInfos ()

11.14.2.4 GetDeviceGroupChannelInfos() [2/2] List<CDeviceGroupChannelInfoTemplateNet^> ^ Get↔
DeviceGroupChannelInfos (
 uint32_t virtualDevice)

11.14.2.5 GetEnabledChannelsInGroup() [1/2] List<bool> ^ GetEnabledChannelsInGroup (
 DacqGroupChannelEnumTemplateNet GroupID)

11.14.2.6 GetEnabledChannelsInGroup() [2/2] List<bool> ^ GetEnabledChannelsInGroup (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

11.14.2.7 GetGroupID() [1/2] DacqGroupChannelEnumTemplateNet GetGroupID (
 uint32_t Index)

11.14.2.8 GetGroupID() [2/2] DacqGroupChannelEnumTemplateNet GetGroupID (
 uint32_t Index,
 uint32_t virtualDevice)

11.14.2.9 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
 DacqGroupChannelEnumTemplateNet GroupID)

11.14.2.10 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

11.14.2.11 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
 DacqGroupChannelEnumTemplateNet GroupID)

11.14.2.12 GetGroupSampleSize() [2/2] [SampleSizeNet](#) GetGroupSampleSize (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

11.14.2.13 GetGroupType() [1/2] [DacqMeaGroupTypeEnumNet](#) GetGroupType (
 DacqGroupChannelEnumTemplateNet GroupID)


11.14.2.14 GetGroupType() [2/2] [DacqMeaGroupTypeEnumNet](#) GetGroupType (
 DacqGroupChannelEnumTemplateNet GroupID,
 uint32_t virtualDevice)

11.14.2.15 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ()

11.14.2.16 GetNumberOfSupportedGroups() [2/2] uint32_t GetNumberOfSupportedGroups (
 uint32_t virtualDevice)

11.15 CDeviceGroupChannelInfoGenericNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoGenericNet:

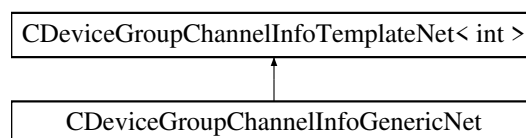

Public Member Functions

- [CDeviceGroupChannelInfoGenericNet](#) (int id, int channels, [DacqMeaGroupTypeEnumNet](#) type)


Additional Inherited Members

11.15.1 Constructor & Destructor Documentation

11.15.1.1 CDeviceGroupChannelInfoGenericNet() `CDeviceGroupChannelInfoGenericNet` (
 `int id,`
 `int channels,`
 `DacqMeaGroupTypeEnumNet type`)

11.16 CDeviceGroupChannelInfoMEA2100_256Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoMEA2100_256Net:

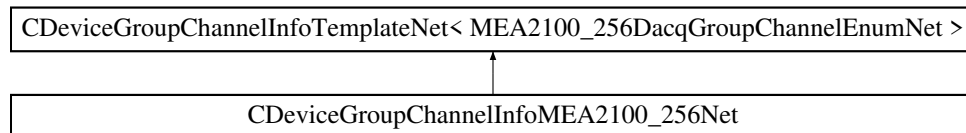

Public Member Functions

- `CDeviceGroupChannelInfoMEA2100_256Net` (`MEA2100_256DacqGroupChannelEnumNet id`, `int channels`, `DacqMeaGroupTypeEnumNet type`)


Additional Inherited Members

11.16.1 Constructor & Destructor Documentation

11.16.1.1 CDeviceGroupChannelInfoMEA2100_256Net() `CDeviceGroupChannelInfoMEA2100_256Net` (
 `MEA2100_256DacqGroupChannelEnumNet id,`
 `int channels,`
 `DacqMeaGroupTypeEnumNet type`)

11.17 CDeviceGroupChannelInfoNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoNet:

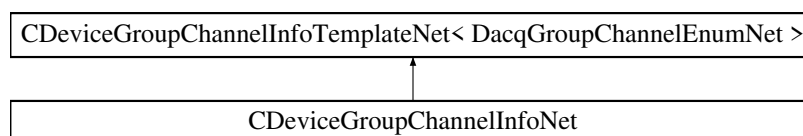

Public Member Functions

- `CDeviceGroupChannelInfoNet` (`DacqGroupChannelEnumNet id`, `int channels`, `DacqMeaGroupTypeEnumNet type`)


Additional Inherited Members

11.17.1 Constructor & Destructor Documentation

11.17.1.1 CDeviceGroupChannelInfoNet() `CDeviceGroupChannelInfoNet` (
 `DacqGroupChannelEnumNet` *id*,
 `int` *channels*,
 `DacqMeaGroupTypeEnumNet` *type*)

11.18 CDeviceGroupChannelInfoSCUNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoSCUNet:

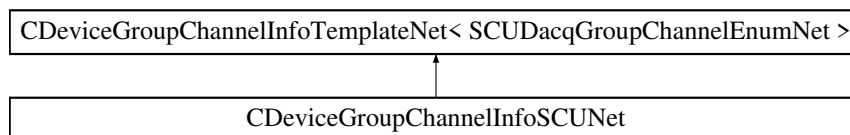

Public Member Functions

- `CDeviceGroupChannelInfoSCUNet` (`SCUDacqGroupChannelEnumNet` *id*, `int` *channels*, `DacqMeaGroupTypeEnumNet` *type*)

Additional Inherited Members

11.18.1 Constructor & Destructor Documentation

11.18.1.1 CDeviceGroupChannelInfoSCUNet() `CDeviceGroupChannelInfoSCUNet` (
 `SCUDacqGroupChannelEnumNet` *id*,
 `int` *channels*,
 `DacqMeaGroupTypeEnumNet` *type*)

11.19 CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet > Class Template Reference

Public Member Functions

- `CDeviceGroupChannelInfoTemplateNet` (`DacqGroupChannelEnumTemplateNet` *id*, `int` *channels*, `DacqMeaGroupTypeEnumNet` *type*)

Public Attributes

- DacqGroupChannelEnumTemplateNet [GroupID](#)
- int [NumberOfChannels](#)
- DacqMeaGroupTypeEnumNet [GroupType](#)

11.19.1 Constructor & Destructor Documentation

11.19.1.1 CDeviceGroupChannelInfoTemplateNet() [CDeviceGroupChannelInfoTemplateNet](#) (
 DacqGroupChannelEnumTemplateNet *id*,
 int *channels*,
 DacqMeaGroupTypeEnumNet *type*)

11.19.2 Member Data Documentation


11.19.2.1 GroupID [DacqGroupChannelEnumTemplateNet](#) [GroupID](#)

11.19.2.2 GroupType [DacqMeaGroupTypeEnumNet](#) [GroupType](#)

11.19.2.3 NumberOfChannels [int](#) [NumberOfChannels](#)

11.20 CDeviceGroupChannelInfoW2100Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoW2100Net:

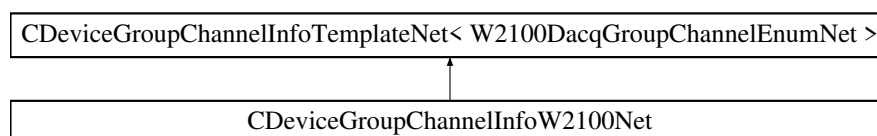**Public Member Functions**

- [CDeviceGroupChannelInfoW2100Net](#) ([W2100DacqGroupChannelEnumNet](#) *id*, int *channels*, [DacqMeaGroupTypeEnumNet](#) *type*)

Additional Inherited Members


11.20.1 Constructor & Destructor Documentation

11.20.1.1 CDeviceGroupChannelInfoW2100Net() `CDeviceGroupChannelInfoW2100Net (W2100DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)`

11.21 CDigOutStimulatorFunctionNet Class Reference

`CDigOutStimulatorFunctionNet` is the class of the DigOut stimulator function class.

Inheritance diagram for `CDigOutStimulatorFunctionNet`:

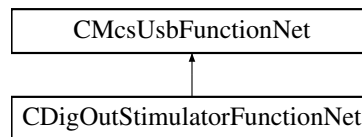

Public Member Functions

- `CDigOutStimulatorFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pDigOutStimulatorFunctionPointerContainer)`
Initializes a new instance of the `CDigOutStimulatorFunctionNet` class.
- `CDigOutStimulatorFunctionNet (CMcsUsbNet^ mcsusb)`
- virtual `~CDigOutStimulatorFunctionNet ()`
- `!CDigOutStimulatorFunctionNet ()`
- void `ClearChannel (int32_t NrChannel)`
clear stimulation pattern
- `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareChannelData (array< int32_t >^ Amplitude, array< uint64_t >^ Duration)`
prepares the channel data for the device and unrolles the data for the GUI
- void `SendChannelData (int32_t NrChannel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled)`
send or append stimulation pattern
- int32_t `GetNumberOfChannels ()`
get the number of channels available on the device
- void `SetGlobalRepeat (int32_t NrChannel, bool value)`
set repeat whole stimulation pattern
- bool `GetGlobalRepeat (int32_t NrChannel)`
get repeat whole stimulation pattern
- void `SetStartTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)`
sets start condition of digital out stimulator
- `DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (int32_t NrChannel)`
queries start condition of digital out stimulator
- void `SetStopTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)`
sets stop condition of digital out stimulator
- `DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope (int32_t NrChannel)`
queries stop condition of digital out stimulator

Additional Inherited Members

11.21.1 Detailed Description

[CDigOutStimulatorFunctionNet](#) is the class of the DigOut stimulator function class.

11.21.2 Constructor & Destructor Documentation

11.21.2.1 CDigOutStimulatorFunctionNet() [1/2] [CDigOutStimulatorFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *pDigOutStimulatorFunctionPointerContainer*)

Initializes a new instance of the [CDigOutStimulatorFunctionNet](#) class.

11.21.2.2 CDigOutStimulatorFunctionNet() [2/2] [CDigOutStimulatorFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

11.21.2.3 ~CDigOutStimulatorFunctionNet() virtual [~CDigOutStimulatorFunctionNet](#) () [virtual]

11.21.2.4 !CDigOutStimulatorFunctionNet() [!CDigOutStimulatorFunctionNet](#) ()

11.21.3 Member Function Documentation

11.21.3.1 ClearChannel() void ClearChannel (
int32_t *NrChannel*)

clear stimulation pattern

Parameters

<i>NrChannel</i>	the channel to clear
------------------	----------------------

11.21.3.2 GetGlobalRepeat() `bool GetGlobalRepeat (`
`int32_t NrChannel)`

get repeat whole stimulation pattern

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

current value

11.21.3.3 GetNumberOfChannels() `int32_t GetNumberOfChannels ()`

get the number of channels available on the device

Returns

the number of channels

11.21.3.4 GetStartTriggerSlope() `DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (`
`int32_t NrChannel)`

queries start condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

start condition (rising or falling edge)

11.21.3.5 GetStopTriggerSlope() `DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope (`
`int32_t NrChannel)`

queries stop condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
------------------	----------------

Returns

stop condition (rising or falling edge)

11.21.3.6 PrepareChannelData() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareChannelData (`
`array< int32_t >^ Amplitude,`
`array< uint64_t >^ Duration)`

prepares the channel data for the device and unrolles the data for the GUI

Parameters

<i>Amplitude</i>	array of amplitudes
<i>Duration</i>	array of durations

Returns

11.21.3.7 SendChannelData() `void SendChannelData (`
`int32_t NrChannel,`
`CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled`
`)`

send or append stimulation pattern

Parameters

<i>NrChannel</i>	the channel to send data to
<i>device_data_and_unrolled</i>	internal, use value obtained from PrepareChannelData

11.21.3.8 SetGlobalRepeat() `void SetGlobalRepeat (`
`int32_t NrChannel,`
`bool value)`

set repeat whole stimulation pattern

Parameters

<i>NrChannel</i>	channel number
<i>value</i>	new value

11.21.3.9 SetStartTriggerSlope() `void SetStartTriggerSlope (`
`int32_t NrChannel,`
`DigitalStimulatorTriggerSlopeEnumNet Condition)`

sets start condition of digital out stimulator

Parameters

<i>NrChannel</i>	channel number
<i>Condition</i>	start condition (rising or falling edge)

11.21.3.10 SetStopTriggerSlope() `void SetStopTriggerSlope (`
`int32_t NrChannel,`
`DigitalStimulatorTriggerSlopeEnumNet Condition)`

sets stop condition of digital out stimulator


Parameters

<i>NrChannel</i>	channel number
<i>Condition</i>	stop condition (rising or falling edge)

11.22 CEncapsulatorDeviceNet Class Reference

[CEncapsulatorDeviceNet](#) is the to control the MCS HiClamp device

Inheritance diagram for CEncapsulatorDeviceNet:

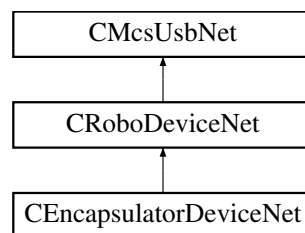

Public Member Functions

- [CEncapsulatorDeviceNet](#) (void)
- [CRoboFluidDeviceNet](#) ^ [GetRoboFluidDevice](#) ()

Additional Inherited Members

11.22.1 Detailed Description

[CEncapsulatorDeviceNet](#) is the to control the MCS HiClamp device

11.22.2 Constructor & Destructor Documentation

11.22.2.1 CEncapsulatorDeviceNet() `CEncapsulatorDeviceNet (void)`


11.22.3 Member Function Documentation

11.22.3.1 GetRoboFluidDevice() `CRoboFluidDeviceNet ^ GetRoboFluidDevice ()`

11.23 CExternDTesterDeviceNet Class Reference

[CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)

Inheritance diagram for CExternDTesterDeviceNet:

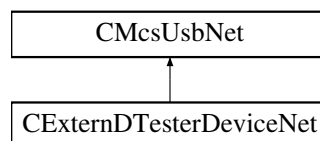

Public Member Functions

- [CExternDTesterDeviceNet](#) ()
Initializes a new instance of the [CExternDTesterDeviceNet](#) class.
- virtual [~CExternDTesterDeviceNet](#) ()
- [!CExternDTesterDeviceNet](#) ()
- `array< uint8_t > ^ Read (int configString_Length)`
Reads the config string from the device.
- `String ^ Read2 ()`
Reads the config string from the device.
- void [Write](#) (array< uint8_t > ^ configString)
Reads the config string from the device.
- void [Write2](#) (String ^ configString)
Reads the config string from the device.

Additional Inherited Members

11.23.1 Detailed Description

[CExternDTesterDeviceNet](#) is the class to access the ExternD Tester (Handheld Device Tester D)

11.23.2 Constructor & Destructor Documentation

11.23.2.1 CExternDTesterDeviceNet() `CExternDTesterDeviceNet ()`

Initializes a new instance of the `CExternDTesterDeviceNet` class.

11.23.2.2 ~CExternDTesterDeviceNet() `virtual ~CExternDTesterDeviceNet () [virtual]`

11.23.2.3 "!CExternDTesterDeviceNet() `!CExternDTesterDeviceNet ()`

11.23.3 Member Function Documentation

11.23.3.1 Read() `array<uint8_t> ^ Read (int configString_Length)`

Reads the config string from the device.

Parameters

<code>configString_Length</code>	The maximal length of configString.
----------------------------------	-------------------------------------

Returns

The config string.

11.23.3.2 Read2() `String ^ Read2 ()`

Reads the config string from the device.

Returns

The config string.

11.23.3.3 Write() `void Write (array< uint8_t > ^ configString)`

Reads the config string from the device.

Parameters

<code>configString</code>	The config string.
---------------------------	--------------------

11.23.3.4 Write2() `void Write2 (String^ configString)`

Reads the config string from the device.

Parameters

<code>configString</code>	The config string.
---------------------------	--------------------

11.24 CFilterCoefficientsNet Class Reference

Classes

- struct [s_FilterAttributesNet](#)

Public Member Functions

- [CFilterCoefficientsNet](#) ()
- [CFilterCoefficientsNet](#) (double b0, double b1, double b2, double a1, double a2)
- [CFilterCoefficientsNet](#) (double b0, double b1, double a1)
- [CFilterCoefficientsNet](#) (array< double >^ b, array< double >^ a)
- [~CFilterCoefficientsNet](#) ()
- bool [IsEqual](#) (array< uint64_t >^ coefficients, [s_FilterAttributesNet](#)^ FiltAttr)
- bool [IsEqual](#) (array< uint64_t >^ coefficients, [s_FilterAttributesNet](#)^ FiltAttr, bool DoMCSLegacyCompare)
- uint64_t [GetUIntB](#) (int index, [s_FilterAttributesNet](#)^ FiltAttr)
- uint64_t [GetUIntA](#) (int index, [s_FilterAttributesNet](#)^ FiltAttr)

Properties

- array< double >^ [A](#) [get]
- array< double >^ [B](#) [get]

11.24.1 Constructor & Destructor Documentation

11.24.1.1 CFilterCoefficientsNet() [1/4] [CFilterCoefficientsNet](#) ()

11.24.1.2 CFilterCoefficientsNet() [2/4] `CFilterCoefficientsNet` (
 double *b0*,
 double *b1*,
 double *b2*,
 double *a1*,
 double *a2*)

11.24.1.3 CFilterCoefficientsNet() [3/4] `CFilterCoefficientsNet` (
 double *b0*,
 double *b1*,
 double *a1*)

11.24.1.4 CFilterCoefficientsNet() [4/4] `CFilterCoefficientsNet` (
 array< double >^ *b*,
 array< double >^ *a*)

11.24.1.5 ~CFilterCoefficientsNet() `~CFilterCoefficientsNet` ()

11.24.2 Member Function Documentation

11.24.2.1 GetUintA() `uint64_t GetUintA` (
 int *index*,
 s_FilterAttributesNet^ *FiltAttr*)

11.24.2.2 GetUintB() `uint64_t GetUintB` (
 int *index*,
 s_FilterAttributesNet^ *FiltAttr*)

11.24.2.3 IsEqual() [1/2] `bool IsEqual` (
 array< uint64_t >^ *coefficients*,
 s_FilterAttributesNet^ *FiltAttr*)

11.24.2.4 IsEqual() [2/2] `bool IsEqual (`
`array< uint64_t >^ coefficients,`
`s_FilterAttributesNet^ FiltAttr,`
`bool DoMCSLegacyCompare)`


11.24.3 Property Documentation

11.24.3.1 A `array< double>^ A [get]`

11.24.3.2 B `array< double>^ B [get]`

11.25 CFilterConfigurationNet Class Reference

Inheritance diagram for CFilterConfigurationNet:

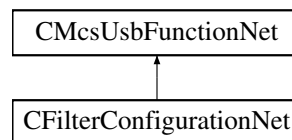

Public Member Functions

- `CFilterConfigurationNet (CMcsUsbNet^ mcsusb)`
- `void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet^ Coefficients, CFilterPropertyNet^ FilterProp)`
- `void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet^ CoefficientsSet1, CFilterCoefficientsNet^ CoefficientsSet2, CFilterPropertyNet^ FilterProp)`
- `void SetFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`
- `void EraseFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`
- `void SetHighpassFilterEnable (bool enable)`
- `bool GetHighpassFilterEnable ()`
- `void ResetHighpassFilter ()`
- `uint32_t GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, FilterAttributeEnumNet index)`
- `CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`

Additional Inherited Members

11.25.1 Constructor & Destructor Documentation

11.25.1.1 CFilterConfigurationNet() `CFilterConfigurationNet (CMcsUsbNet^ mcsusb)`

11.25.2 Member Function Documentation

11.25.2.1 EraseFilterParameterPermanent() `void EraseFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`

11.25.2.2 GetFilterAttributes() [1/2] `CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`

11.25.2.3 GetFilterAttributes() [2/2] `uint32_t GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, FilterAttributeEnumNet index)`

11.25.2.4 GetHighpassFilterEnable() `bool GetHighpassFilterEnable ()`

11.25.2.5 ResetHighpassFilter() `void ResetHighpassFilter ()`

11.25.2.6 SetFilterParameter() [1/2] `void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet^ Coefficients, CFilterPropertyNet^ FilterProp)`


11.25.2.7 SetFilterParameter() [2/2] `void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet^ CoefficientsSet1, CFilterCoefficientsNet^ CoefficientsSet2, CFilterPropertyNet^ FilterProp)`

11.25.2.8 SetFilterParameterPermanent() `void SetFilterParameterPermanent (`
`DacqGroupChannelEnumNet GroupID,`
`uint32_t FilterNumber)`

11.25.2.9 SetHighpassFilterEnable() `void SetHighpassFilterEnable (`
`bool enable)`

11.26 CFilterConfigurationRegisterNet Class Reference

Inheritance diagram for CFilterConfigurationRegisterNet:

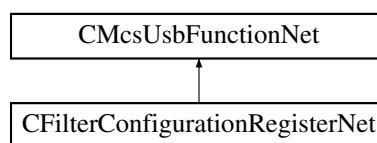

Public Member Functions

- [CFilterConfigurationRegisterNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetFilterParameter](#) (uint32_t FilterCoefRegBase, [CFilterCoefficientsNet](#)[^] Coefficients, uint32_t FilterInfoRegBase, [CFilterPropertyNet](#)[^] FilterProp)
- void [SetFilterParameter](#) (uint32_t FilterCoefSet1RegBase, [CFilterCoefficientsNet](#)[^] CoefficientsSet1, uint32_t FilterCoefSet2RegBase, [CFilterCoefficientsNet](#)[^] CoefficientsSet2, uint32_t FilterInfoRegBase, [CFilterPropertyNet](#)[^] FilterProp)
- void [SetFilterParameterPermanent](#) (uint32_t FilterCoefRegBase, uint32_t FilterCoefDmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void [SetFilterParameterPermanent](#) (uint32_t FilterCoefSet1RegBase, uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2RegBase, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void [EraseFilterParameterPermanent](#) (uint32_t FilterCoefDmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void [EraseFilterParameterPermanent](#) (uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)

Additional Inherited Members

11.26.1 Constructor & Destructor Documentation

11.26.1.1 CFilterConfigurationRegisterNet() `CFilterConfigurationRegisterNet (`
`CMcsUsbNet^ mcsusb)`

11.26.2 Member Function Documentation

11.26.2.1 EraseFilterParameterPermanent() [1/2] `void EraseFilterParameterPermanent (`
 `uint32_t FilterCoefDmaReg,`
 `uint32_t FilterInfoDmaReg,`
 `uint32_t EEPROMBase,`
 `uint32_t EEPROMSize)`

11.26.2.2 EraseFilterParameterPermanent() [2/2] `void EraseFilterParameterPermanent (`
 `uint32_t FilterCoefSet1DmaReg,`
 `uint32_t FilterCoefSet2DmaReg,`
 `uint32_t FilterInfoDmaReg,`
 `uint32_t EEPROMBase,`
 `uint32_t EEPROMSize)`

11.26.2.3 SetFilterParameter() [1/2] `void SetFilterParameter (`
 `uint32_t FilterCoefRegBase,`
 `CFilterCoefficientsNet^ Coefficients,`
 `uint32_t FilterInfoRegBase,`
 `CFilterPropertyNet^ FilterProp)`

11.26.2.4 SetFilterParameter() [2/2] `void SetFilterParameter (`
 `uint32_t FilterCoefSet1RegBase,`
 `CFilterCoefficientsNet^ CoefficientsSet1,`
 `uint32_t FilterCoefSet2RegBase,`
 `CFilterCoefficientsNet^ CoefficientsSet2,`
 `uint32_t FilterInfoRegBase,`
 `CFilterPropertyNet^ FilterProp)`

11.26.2.5 SetFilterParameterPermanent() [1/2] `void SetFilterParameterPermanent (`
 `uint32_t FilterCoefRegBase,`
 `uint32_t FilterCoefDmaReg,`
 `uint32_t FilterInfoRegBase,`
 `uint32_t FilterInfoDmaReg,`
 `uint32_t EEPROMBase,`
 `uint32_t EEPROMSize)`

11.26.2.6 SetFilterParameterPermanent() [2/2] `void SetFilterParameterPermanent (`
 `uint32_t FilterCoefSet1RegBase,`
 `uint32_t FilterCoefSet1DmaReg,`
 `uint32_t FilterCoefSet2RegBase,`
 `uint32_t FilterCoefSet2DmaReg,`
 `uint32_t FilterInfoRegBase,`
 `uint32_t FilterInfoDmaReg,`
 `uint32_t EEPROMBase,`
 `uint32_t EEPROMSize)`

11.27 CFilterPropertyNet Class Reference

Public Member Functions

- [CFilterPropertyNet](#) (uint32_t *CornerFrequenzymHz*, uint32_t *Order*, [FilterBandEnumNet](#) *FilterBand*, [FilterFamilyEnumNet](#) *FilterFamily*, [FilterTypeEnumNet](#) *FilterType*, bool *Active*)
- [~CFilterPropertyNet](#) ()
- virtual `System::String ^ ToString ()` override

Properties

- uint32_t [CornerFrequencymHz](#) [get]
- double [CornerFrequency](#) [get]
- uint32_t [Order](#) [get]
- [FilterBandEnumNet](#) *FilterBand* [get]
- [FilterFamilyEnumNet](#) *FilterFamily* [get]
- [FilterTypeEnumNet](#) *FilterType* [get]
- bool [FilterActive](#) [get]

11.27.1 Constructor & Destructor Documentation

11.27.1.1 CFilterPropertyNet() [CFilterPropertyNet](#) (
 uint32_t *CornerFrequenzymHz*,
 uint32_t *Order*,
 [FilterBandEnumNet](#) *FilterBand*,
 [FilterFamilyEnumNet](#) *FilterFamily*,
 [FilterTypeEnumNet](#) *FilterType*,
 bool *Active*)

11.27.1.2 ~CFilterPropertyNet() [~CFilterPropertyNet](#) ()

11.27.2 Member Function Documentation

11.27.2.1 ToString() virtual `System::String ^ ToString ()` [override], [virtual]

11.27.3 Property Documentation

11.27.3.1 CornerFrequency `double CornerFrequency [get]`

11.27.3.2 CornerFrequencymHz `uint32_t CornerFrequencymHz [get]`

11.27.3.3 FilterActive `bool FilterActive [get]`

11.27.3.4 FilterBand `FilterBandEnumNet FilterBand [get]`

11.27.3.5 FilterFamily `FilterFamilyEnumNet FilterFamily [get]`


11.27.3.6 FilterType `FilterTypeEnumNet FilterType [get]`

11.27.3.7 Order `uint32_t Order [get]`

11.28 CFluidControlDeviceNet Class Reference

[CFluidControlDeviceNet](#) is the class to control MCS FluidControl (FCB and FCX) device.

Inheritance diagram for CFluidControlDeviceNet:

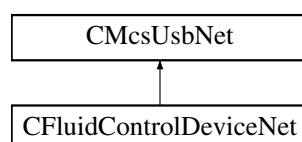

Public Member Functions

- [CFluidControlDeviceNet](#) ()
Initialize a new instance of the [CFluidControlDeviceNet](#) class.
- [~CFluidControlDeviceNet](#) ()
Default destructor.
- void [SetValve](#) (unsigned int value)
Open or Close valves.
- void [SetSingleValve](#) (unsigned short valve, unsigned short onoff)
Opens or Closes a valve.
- void [SetDigout](#) (unsigned int value)
Define the pattern on the Digital Output.
- void [SetPWM](#) (unsigned int channel, unsigned int value)
Sets the duty cycle of the PWM output.
- void [CalibrateThermocouple](#) (unsigned int channel)
Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.
- void [SetThermocoupleNanovoltPerKelvin](#) (unsigned int channel, unsigned int value)
Sets the proportional constant for the Thermocouple.
- unsigned int [GetValve](#) ()
Gets the state of the valves.
- unsigned short [GetSingleValve](#) (unsigned short valve)
Gets the state of a valve.
- unsigned int [GetDigout](#) ()
Gets the state of the digital output.
- unsigned int [GetPWM](#) (unsigned int channel)
Gets the state of the PWM output.
- unsigned int [GetAdc](#) (unsigned int channel)
Reads an ADC Value.
- unsigned int [GetDigin](#) ()
Reads the digital input.
- int [GetThermocoupleTemperature](#) (unsigned int channel)
Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.
- int [GetReferenceTemperature](#) (unsigned int channel)
Reads the reference temperature for the Thermocouple.
- unsigned int [GetThermocoupleCalibration](#) (unsigned int channel)
Gets the calibration constant for the Thermocouple ADC.
- unsigned int [GetThermocoupleNanovoltPerKelvin](#) (unsigned int channel)
Reads the proportional constant for the Thermocouple.

Properties

- [CMcsBus_VoltageModeNet](#)[^] [McsBus_VoltageMode](#) [get]

Additional Inherited Members

11.28.1 Detailed Description

[CFluidControlDeviceNet](#) is the class to control MCS FluidControl (FCB and FCX) device.

11.28.2 Constructor & Destructor Documentation

11.28.2.1 CFluidControlDeviceNet() `CFluidControlDeviceNet ()`

Initialize a new instance of the `CFluidControlDeviceNet` class.

11.28.2.2 ~CFluidControlDeviceNet() `~CFluidControlDeviceNet ()`

Default destructor.

11.28.3 Member Function Documentation

11.28.3.1 CalibrateThermocouple() `void CalibrateThermocouple (unsigned int channel)`

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

11.28.3.2 GetAdc() `unsigned int GetAdc (unsigned int channel)`

Reads an ADC Value.

Parameters

<i>channel</i>	The ADC channel number to query.
----------------	----------------------------------

Returns

The current ADC value.

11.28.3.3 GetDigin() `unsigned int GetDigin ()`

Reads the digital input.

Returns

The bit pattern of the state of the digital inputs.

11.28.3.4 GetDigout() `unsigned int GetDigout ()`

Gets the state of the digital output.

Returns

The current state of the digital outputs as a bit pattern.

11.28.3.5 GetPWM() `unsigned int GetPWM (`
`unsigned int channel)`

Gets the state of the PWM output.

Returns

The current state of the PWM outputs duty cycle in permille.

11.28.3.6 GetReferenceTemperature() `int GetReferenceTemperature (`
`unsigned int channel)`

Reads the reference temperature for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The temperature from the Thermocouple in 1/100 °C.

11.28.3.7 GetSingleValve() `unsigned short GetSingleValve (`
`unsigned short valve)`

Gets the state of a valve.

Parameters

<i>valve</i>	number of valve
--------------	-----------------

Returns

state of the valve

11.28.3.8 GetThermocoupleCalibration() `unsigned int GetThermocoupleCalibration (unsigned int channel)`

Gets the calibration constant for the Thermocouple ADC.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The calibration constant for the Thermocouple ADC.

11.28.3.9 GetThermocoupleNanovoltPerKelvin() `unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)`

Reads the proportional constant for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

11.28.3.10 GetThermocoupleTemperature() `int GetThermocoupleTemperature (unsigned int channel)`

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The temperature difference between both Thermocouple junctions in 1/100 °C.

11.28.3.11 GetValve() `unsigned int GetValve ()`

Gets the state of the valves.

Returns

The current state of the valves as a bit pattern.

11.28.3.12 SetDigout() `void SetDigout (`
`unsigned int value)`

Define the pattern on the Digital Output.

Parameters

<i>value</i>	bit pattern on the digital output.
--------------	------------------------------------

11.28.3.13 SetPWM() `void SetPWM (`
`unsigned int channel,`
`unsigned int value)`

Sets the duty cycle of the PWM output.

Parameters

<i>channel</i>	PWM channel number.
<i>value</i>	duty cycle of the PWM output in permille.

11.28.3.14 SetSingleValve() `void SetSingleValve (`
`unsigned short valve,`
`unsigned short onoff)`

Opens or Closes a valve.

Parameters

<i>valve</i>	number of valve to be changed.
--------------	--------------------------------

Parameters

<i>onoff</i>	open or close the valve.
--------------	--------------------------

11.28.3.15 SetThermocoupleNanovoltPerKelvin() `void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)`

Sets the proportinal constant for the Thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
<i>value</i>	proportional constant for the Thermocouple in Nanovolt per Kelvin.

11.28.3.16 SetValve() `void SetValve (unsigned int value)`

Open or Close valves.

Parameters

<i>value</i>	bit pattern of valves which should be open.
--------------	---


11.28.4 Property Documentation

11.28.4.1 McsBus_VoltageMode `CMcsBus_VoltageModeNet^ McsBus_VoltageMode [get]`

11.29 CFYIDeviceNet Class Reference

[CFYIDeviceNet](#) is the class to control the MCS FYI device

Inheritance diagram for CFYIDeviceNet:

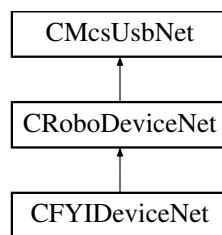

Public Member Functions

- [CFYIDeviceNet](#) (void)

Properties

- [CRobo_FYITemp_FunctionNet](#)^ [FYITemp](#) [get]
- [CRobo_FYIProgram_FunctionNet](#)^ [FYIProgram](#) [get]
- [CMcsBus_SensorNet](#)^ [Sensor](#) [get]

Additional Inherited Members

11.29.1 Detailed Description

[CFYIDeviceNet](#) is the class to control the MCS FYI device

11.29.2 Constructor & Destructor Documentation

11.29.2.1 [CFYIDeviceNet\(\)](#) [CFYIDeviceNet](#) (
void)

11.29.3 Property Documentation

11.29.3.1 [FYIProgram](#) [CRobo_FYIProgram_FunctionNet](#)^ [FYIProgram](#) [get]


11.29.3.2 [FYITemp](#) [CRobo_FYITemp_FunctionNet](#)^ [FYITemp](#) [get]

11.29.3.3 Sensor `CMcsBus_SensorNet`^ Sensor [get]

11.30 CGenericDevelopDeviceNet Class Reference

`CGenericDevelopDeviceNet` is the class to use during development of a new device.

Inheritance diagram for `CGenericDevelopDeviceNet`:

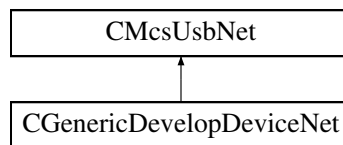

Public Member Functions

- `CGenericDevelopDeviceNet` (void)
Initialize a new instance of the `CGenericDevelopDeviceNet` class.
- `~CGenericDevelopDeviceNet` (void)
- void `SetValue` (uint16_t value, uint16_t index)
Sets .

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

- template<typename C >
void `SetBuffer` (uint16_t value, uint16_t index, array< C >^ buffer)
- void `SetUByteBuffer` (uint16_t value, uint16_t index, array< unsigned char >^ buffer)
Sends an array of type unsigned char to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetByteBuffer](#) (uint16_t value, uint16_t Index, array< char >^ buffer)

Sends an array of type char to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

Index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetUShortBuffer](#) (uint16_t value, uint16_t index, array< unsigned short >^ buffer)

Sends an array of type unsigned short to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetShortBuffer](#) (uint16_t value, uint16_t index, array< short >^ buffer)

Sends an array of type short to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetUIntBuffer](#) (uint16_t value, uint16_t index, array< unsigned int >^ buffer)
Sends an array of unsigned int to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- void [SetIntBuffer](#) (uint16_t value, uint16_t index, array< int >^ buffer)
Sends an array of type int to the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

buffer	<i>The buffer to send.</i>
--------	----------------------------

- `template<typename C >`
`array< C > ^ GetBuffer (uint16_t value, uint16_t index, int size)`
- `array< unsigned char > ^ GetUByteBuffer (uint16_t value, uint16_t index, int size)`
Gets an array of type unsigned char from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- `array< char > ^ GetByteBuffer (uint16_t value, uint16_t index, int size)`
Gets an array of type char from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- `array< unsigned short > ^ GetUShortBuffer (uint16_t value, uint16_t index, int size)`
Gets an array of type unsigned short from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< short > ^ [GetShortBuffer](#) (uint16_t value, uint16_t index, int size)

Gets an array of type short from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- array< unsigned int > ^ [GetUIntBuffer](#) (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned int from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- `array< int > ^ GetIntBuffer (uint16_t value, uint16_t index, int size)`

Gets an array of type int from the device.

Parameters

value	<i>The value of the request.</i>
-------	----------------------------------

Parameters

index	<i>The index of the request.</i>
-------	----------------------------------

Parameters

size	<i>The size of the array.</i>
------	-------------------------------

Returns

The array of data from the device.

- `template<typename C >`
`void VendorOutRequest (uint8_t request, uint16_t value, uint16_t index, array< C >^ buffer)`
- `template<typename C >`
`array< C > ^ VendorInRequest (uint8_t request, uint16_t value, uint16_t index, int size)`
- `array< uint8_t > ^ FindEndpoints (uint8_t type, uint8_t direction)`
- `IntPtr OpenPipe (uint8_t endpointAddress)`

Open a Pipe to an USB Endpoint.

Parameters

endpointAddress	<i>The Endpoint Number.</i>
-----------------	-----------------------------

*Returns**A handle to the endpoint.*

- void [ClosePipe](#) (IntPtr pipeHandle)
Close a Pipe to an USB Endpoint.

Parameters

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

- void [ResetPipe](#) (IntPtr pipeHandle)
Reset a Pipe to an USB Endpoint.

Parameters

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

- template<typename C >
array< C > ^ [ReadPipe](#) (IntPtr pipeHandle, uint32_t size)
Read data from an USB Endpoint.

Parameters

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

Parameters

size	<i>Number of elements to read.</i>
------	------------------------------------

*Returns**An array of data read.*

- template<typename C >
void [WritePipe](#) (IntPtr pipeHandle, array< C >^ buffer)
Write data to an USB Endpoint.

Parameters

pipeHandle	<i>The endpoint handle.</i>
------------	-----------------------------

Parameters

buffer	<i>An array of data to write.</i>
--------	-----------------------------------

Additional Inherited Members

11.30.1 Detailed Description

[CGenericDevelopDeviceNet](#) is the class to use during development of a new device.

11.30.2 Constructor & Destructor Documentation

11.30.2.1 CGenericDevelopDeviceNet() [CGenericDevelopDeviceNet](#) (
void)

Initialize a new instance of the [CGenericDevelopDeviceNet](#) class.

11.30.2.2 ~CGenericDevelopDeviceNet() [~CGenericDevelopDeviceNet](#) (
void)

11.30.3 Member Function Documentation

11.30.3.1 ClosePipe() void ClosePipe (
IntPtr *pipeHandle*)

Close a Pipe to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

11.30.3.2 FindEndpoints() array<uint8_t> ^ FindEndpoints (
uint8_t *type*,
uint8_t *direction*)

11.30.3.3 GetBuffer() `array<C> ^ GetBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

11.30.3.4 GetByteBuffer() `array<char> ^ GetByteBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type char from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

11.30.3.5 GetIntBuffer() `array<int> ^ GetIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type int from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

11.30.3.6 GetShortBuffer() `array<short> ^ GetShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type short from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

11.30.3.7 GetUByteBuffer() `array<unsigned char> ^ GetUByteBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type unsigned char from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

11.30.3.8 GetUIntBuffer() `array<unsigned int> ^ GetUIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

Gets an array of type unsigned int from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

```
11.30.3.9 GetUShortBuffer() array<unsigned short> ^ GetUShortBuffer (
    uint16_t value,
    uint16_t index,
    int size )
```

Gets an array of type unsigned short from the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>size</i>	The size of the array.
-------------	------------------------

Returns

The array of data from the device.

11.30.3.10 OpenPipe() `IntPtr OpenPipe (`
`uint8_t endpointAddress)`

Open a Pipe to an USB Endpoint.

Parameters

<i>endpointAddress</i>	The Endpoint Number.
------------------------	----------------------

Returns

A handle to the endpoint.

11.30.3.11 ReadPipe() `array<C> ^ ReadPipe (`
`IntPtr pipeHandle,`
`uint32_t size)`

Read data from an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

Parameters

<i>size</i>	Number of elements to read.
-------------	-----------------------------

Returns

An array of data read.

11.30.3.12 ResetPipe() `void ResetPipe (`
 `IntPtr pipeHandle)`

Reset a Pipe to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------

11.30.3.13 SetBuffer() `void SetBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< C >^ buffer)`

11.30.3.14 SetByteBuffer() `void SetByteBuffer (`
 `uint16_t value,`
 `uint16_t Index,`
 `array< char >^ buffer)`

Sends an array of type char to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>Index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.15 SetIntBuffer() `void SetIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< int >^ buffer)`

Sends an array of type int to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.16 SetShortBuffer() `void SetShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< short >^ buffer)`

Sends an array of type short to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.17 SetUByteBuffer() `void SetUByteBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< unsigned char >^ buffer)`

Sends an array of type unsigned char to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.18 SetUIntBuffer() `void SetUIntBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< unsigned int >^ buffer)`

Sends an array of unsigned int to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.19 SetUShortBuffer() `void SetUShortBuffer (`
 `uint16_t value,`
 `uint16_t index,`
 `array< unsigned short >^ buffer)`

Sends an array of type unsigned short to the device.

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

Parameters

<i>buffer</i>	The buffer to send.
---------------	---------------------

11.30.3.20 SetValue() `void SetValue (`
 `uint16_t value,`
 `uint16_t index)`

Sets .

Parameters

<i>value</i>	The value of the request.
--------------	---------------------------

Parameters

<i>index</i>	The index of the request.
--------------	---------------------------

11.30.3.21 VendorInRequest() `array<C> ^ VendorInRequest (`
 `uint8_t request,`
 `uint16_t value,`
 `uint16_t index,`
 `int size)`

11.30.3.22 VendorOutRequest() `void VendorOutRequest (`
 `uint8_t request,`
 `uint16_t value,`
 `uint16_t index,`
 `array< C >^ buffer)`

11.30.3.23 WritePipe() `void WritePipe (`
 `IntPtr pipeHandle,`
 `array< C >^ buffer)`

Write data to an USB Endpoint.

Parameters

<i>pipeHandle</i>	The endpoint handle.
-------------------	----------------------


Parameters

<i>buffer</i>	An array of data to write.
---------------	----------------------------

11.31 CGilsonDeviceNet Class Reference

[CGilsonDeviceNet](#) is the class to control a Gilson device.

Inheritance diagram for CGilsonDeviceNet:

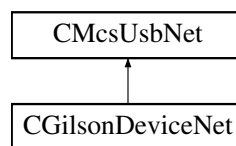

Public Member Functions

- [CGilsonDeviceNet](#) (void)
Initialize a new instance of the [CGilsonDeviceNet](#) class.
- [~CGilsonDeviceNet](#) (void)
- void [ConnectSlave](#) (byte ID)
- void [SendImmediate](#) (wchar_t command)
- String ^ [SendImmediateGetResponse](#) (wchar_t command)
- void [SendBuffered](#) (String^ command)
- String ^ [GetLastAnswer](#) ()

Protected Attributes

- CGilsonDevice * [m_pGilsonDevice](#)

Additional Inherited Members

11.31.1 Detailed Description

[CGilsonDeviceNet](#) is the class to control a Gilson device.

11.31.2 Constructor & Destructor Documentation

11.31.2.1 CGilsonDeviceNet() `CGilsonDeviceNet (`
`void)`

Initialize a new instance of the `CGilsonDeviceNet` class.

11.31.2.2 ~CGilsonDeviceNet() `~CGilsonDeviceNet (`
`void)`

11.31.3 Member Function Documentation

11.31.3.1 ConnectSlave() `void ConnectSlave (`
`byte ID)`

11.31.3.2 GetLastAnswer() `String ^ GetLastAnswer ()`

11.31.3.3 SendBuffered() `void SendBuffered (`
`String^ command)`

11.31.3.4 SendImmediate() `void SendImmediate (`
`wchar_t command)`


11.31.3.5 SendImmediateGetResponse() `String ^ SendImmediateGetResponse (`
`wchar_t command)`

11.31.4 Member Data Documentation

11.31.4.1 m_pGilsonDevice `CGilsonDevice* m_pGilsonDevice [protected]`

11.32 CGrapheneASICDeviceNet Class Reference

Inheritance diagram for CGrapheneASICDeviceNet:

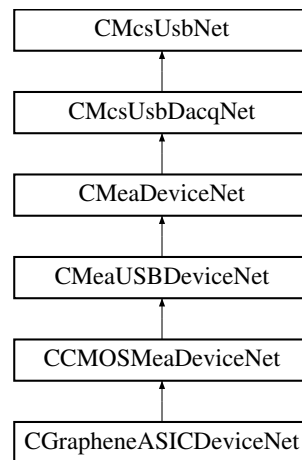

Public Member Functions

- [CGrapheneASICDeviceNet](#) (void)
- [~CGrapheneASICDeviceNet](#) ()
- void [SetBaseSamplerate](#) (int BaseSamplerate) override
- array< int > ^ [GetAvailableBaseSamplerates](#) () override
- void [SetRegionOfInterests](#) (System::Collections::Generic::Dictionary< int, [CCMOSMeaDeviceNet::CRegionOfInterestRect](#)>^ rois) override

Additional Inherited Members

11.32.1 Constructor & Destructor Documentation

11.32.1.1 CGrapheneASICDeviceNet() [CGrapheneASICDeviceNet](#) (void)

11.32.1.2 ~CGrapheneASICDeviceNet() [~CGrapheneASICDeviceNet](#) ()

11.32.2 Member Function Documentation

11.32.2.1 GetAvailableBaseSamplerates() `array<int> ^ GetAvailableBaseSamplerates () [override], [virtual]`

Reimplemented from [CCMOSMeaDeviceNet](#).

11.32.2.2 SetBaseSamplerate() `void SetBaseSamplerate (int BaseSamplerate) [override], [virtual]`


Reimplemented from [CCMOSMeaDeviceNet](#).

11.32.2.3 SetRegionOfInterests() `void SetRegionOfInterests (System::Collections::Generic::Dictionary< int, CCMOSMeaDeviceNet::CRegionOfInterestRect ^>^ rois) [override]`

11.33 CGrapheneFunctionNet Class Reference

[CGrapheneFunctionNet](#) is the class to control Graphene device functions

Inheritance diagram for CGrapheneFunctionNet:

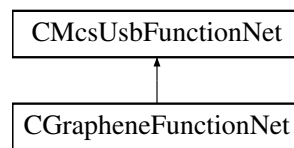

Public Member Functions

- [CGrapheneFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pGrapheneFunctionPointerContainer)
Initializes a new instance of the [CGrapheneFunctionNet](#) class.
- [CGrapheneFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual [~CGrapheneFunctionNet](#) ()
- [ICGrapheneFunctionNet](#) ()
- void [GetVdVsDAC](#) ([System::Runtime::InteropServices::Out]int16_t% Vd, [System::Runtime::InteropServices::Out]int16_t% Vs)
Gets Vd and Vs
- void [GetVdVsDAC](#) (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% Vd, [System::Runtime::InteropServices::Out]int16_t% Vs)
Gets Vd and Vs
- void [SetVdVsDAC](#) (int16_t Vd, int16_t Vs)
Sets Vd and Vs
- void [SetVdVsDAC](#) (uint32_t Headstage, int16_t Vd, int16_t Vs)
Sets Vd and VS
- bool [GetVoltageReached](#) ()
Gets the reached voltage

- bool [GetVoltageReached](#) (uint32_t Headstage)
Gets the reached voltage
- int32_t [GetVoltageRange](#) ()
Gets the voltage range
- int32_t [GetVoltageRange](#) (uint32_t Headstage)
Gets the voltage range
- void [SetVoltageRange](#) (int32_t range)
Sets the voltage range
- void [SetVoltageRange](#) (uint32_t Headstage, int32_t range)
Sets the voltage range
- int32_t [GetVoltageResolution](#) ()
Gets the voltage resolution
- int32_t [GetVoltageResolution](#) (uint32_t Headstage)
Gets the voltage resolution
- void [SetVoltageResolution](#) (int32_t resolution)
Sets the voltage resolution
- void [SetVoltageResolution](#) (uint32_t Headstage, int32_t resolution)
Sets the voltage resolution
- void [GetDACOffset](#) ([System::Runtime::InteropServices::Out]int16_t% offset_vd, [System::Runtime::InteropServices::Out]int16_t% offset_vs)
Gets the DAC offset
- void [GetDACOffset](#) (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% offset_vd, [System::Runtime::InteropServices::Out]int16_t% offset_vs)
Gets the DAC offset
- void [SetDACOffset](#) (int16_t offset_vd, int16_t offset_vs)
Sets the DAC offset
- void [SetDACOffset](#) (uint32_t Headstage, int16_t offset_vd, int16_t offset_vs)
Set the DAC offset
- void [GetVdVs](#) ([System::Runtime::InteropServices::Out]int32_t% Vd, [System::Runtime::InteropServices::Out]int32_t% Vs)
Gets Vd and Vs
- void [GetVdVs](#) (uint32_t Headstage, [System::Runtime::InteropServices::Out]int32_t% Vd, [System::Runtime::InteropServices::Out]int32_t% Vs)
Gets Vd and Vs
- void [SetVdVs](#) (int32_t Vd, int32_t Vs)
Sets Vd and Vs
- void [SetVdVs](#) (uint32_t Headstage, int32_t Vd, int32_t Vs)
Sets Vd and Vs
- void [SetVdsVgs](#) (int32_t Vds, int32_t Vgs)
Sets Vds and Vgs
- void [SetVdsVgs](#) (uint32_t Headstage, int32_t Vds, int32_t Vgs)
Sets Vds and Vgs
- int32_t [GetCur2VolResistance](#) ()
Gets the resistance of the current to voltage converter
- int32_t [GetCur2VolResistance](#) (uint32_t Headstage)
Gets the resistance of the current to voltage converter
- void [GetVdsVgs](#) ([System::Runtime::InteropServices::Out]int32_t% Vds, [System::Runtime::InteropServices::Out]int32_t% Vgs)
Gets Vds and Vgs
- void [GetVdsVgs](#) (uint32_t Headstage, [System::Runtime::InteropServices::Out]int32_t% Vds, [System::Runtime::InteropServices::Out]int32_t% Vgs)
Gets Vds and Vgs

Additional Inherited Members

11.33.1 Detailed Description

[CGrapheneFunctionNet](#) is the class to control Graphene device functions

11.33.2 Constructor & Destructor Documentation

11.33.2.1 CGrapheneFunctionNet() [1/2] [CGrapheneFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pGrapheneFunctionPointerContainer*)

Initializes a new instance of the [CGrapheneFunctionNet](#) class.

11.33.2.2 CGrapheneFunctionNet() [2/2] [CGrapheneFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

11.33.2.3 ~CGrapheneFunctionNet() [virtual ~CGrapheneFunctionNet](#) () [virtual]

11.33.2.4 !CGrapheneFunctionNet() [!CGrapheneFunctionNet](#) ()

11.33.3 Member Function Documentation

11.33.3.1 GetCur2VolResistance() [1/2] [int32_t GetCur2VolResistance](#) ()

Gets the resistance of the current to voltage converter

Returns

The resistance in Ohm

11.33.3.2 GetCur2VolResistance() [2/2] [int32_t GetCur2VolResistance](#) (
 [uint32_t Headstage](#))

Gets the resistance of the current to voltage converter

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The resistance in Ohm

11.33.3.3 GetDACOffset() [1/2] void GetDACOffset (
[System::Runtime::InteropServices::Out] int16_t% *offset_vd*,
[System::Runtime::InteropServices::Out] int16_t% *offset_vs*)

Gets the DAC offset

Parameters

<i>offset_vd</i>	Vd offset in DAC Units
<i>offset_vs</i>	Vs offset in DAC Units

11.33.3.4 GetDACOffset() [2/2] void GetDACOffset (
uint32_t *Headstage*,
[System::Runtime::InteropServices::Out] int16_t% *offset_vd*,
[System::Runtime::InteropServices::Out] int16_t% *offset_vs*)

Gets the DAC offset

Parameters

<i>Headstage</i>	The headstage to query.
<i>offset_vd</i>	Vd offset in DAC Units
<i>offset_vs</i>	Vs offset in DAC Units

11.33.3.5 GetVdsVgs() [1/2] void GetVdsVgs (
[System::Runtime::InteropServices::Out] int32_t% *Vds*,
[System::Runtime::InteropServices::Out] int32_t% *Vgs*)

Gets Vds and Vgs

Parameters

<i>Vds</i>	Vds in μ V
<i>Vgs</i>	Vgs in μ V

11.33.3.6 GetVdsVgs() [2/2] void GetVdsVgs (
 uint32_t *Headstage*,
 [System::Runtime::InteropServices::Out] int32_t% *Vds*,
 [System::Runtime::InteropServices::Out] int32_t% *Vgs*)

Gets Vds and Vgs

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vds</i>	Vds in μV
<i>Vgs</i>	Vgs in μV

11.33.3.7 GetVdVs() [1/2] void GetVdVs (
 [System::Runtime::InteropServices::Out] int32_t% *Vd*,
 [System::Runtime::InteropServices::Out] int32_t% *Vs*)

Gets Vd and Vs

Parameters

<i>Vd</i>	Vd in μV
<i>Vs</i>	Vs in μV

11.33.3.8 GetVdVs() [2/2] void GetVdVs (
 uint32_t *Headstage*,
 [System::Runtime::InteropServices::Out] int32_t% *Vd*,
 [System::Runtime::InteropServices::Out] int32_t% *Vs*)

Gets Vd and Vs

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vd</i>	Vd in μV
<i>Vs</i>	Vs in μV

11.33.3.9 GetVdVsDAC() [1/2] void GetVdVsDAC (
 [System::Runtime::InteropServices::Out] int16_t% *Vd*,
 [System::Runtime::InteropServices::Out] int16_t% *Vs*)

Gets Vd and Vs

Parameters

<i>Vd</i>	Vd in DAC Units
<i>Vs</i>	Vs in DAC Units

11.33.3.10 GetVdVsDAC() [2/2] void GetVdVsDAC (
 uint32_t *Headstage*,
 [System::Runtime::InteropServices::Out] int16_t% *Vd*,
 [System::Runtime::InteropServices::Out] int16_t% *Vs*)

Gets Vd and Vs

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vd</i>	Vd in DAC Units
<i>Vs</i>	Vs in DAC Units

11.33.3.11 GetVoltageRange() [1/2] int32_t GetVoltageRange ()

Gets the voltage range

Returns

The voltage range in mV

11.33.3.12 GetVoltageRange() [2/2] int32_t GetVoltageRange (
 uint32_t *Headstage*)

Gets the voltage range

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The voltage range in mV

11.33.3.13 GetVoltageReached() [1/2] `bool GetVoltageReached ()`

Gets the reached voltage

Returns

the reached voltage

11.33.3.14 GetVoltageReached() [2/2] `bool GetVoltageReached (
uint32_t Headstage)`

Gets the reached voltage

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The reached voltage

11.33.3.15 GetVoltageResolution() [1/2] `int32_t GetVoltageResolution ()`

Gets the voltage resolution

Returns

The voltage resolution in $\mu\text{V}/\text{digit}$

11.33.3.16 GetVoltageResolution() [2/2] `int32_t GetVoltageResolution (
uint32_t Headstage)`

Gets the voltage resolution

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The voltage resolution in $\mu\text{V}/\text{digit}$

11.33.3.17 SetDACOffset() [1/2] `void SetDACOffset (`
 `int16_t offset_vd,`
 `int16_t offset_vs)`

Sets the DAC offset

Parameters

<i>offset_vd</i>	Vd offset in DAC Units
<i>offset_vs</i>	Vs offset in DAC Units

11.33.3.18 SetDACOffset() [2/2] `void SetDACOffset (`
 `uint32_t Headstage,`
 `int16_t offset_vd,`
 `int16_t offset_vs)`

Set the DAC offset

Parameters

<i>Headstage</i>	The headstage to query.
<i>offset_vd</i>	Vd offset in DAC Units
<i>offset_vs</i>	Vs offset in DAC Units

11.33.3.19 SetVdsVgs() [1/2] `void SetVdsVgs (`
 `int32_t Vds,`
 `int32_t Vgs)`

Sets Vds and Vgs

Parameters

<i>Vds</i>	Vds in μV
<i>Vgs</i>	Vgs in μV

11.33.3.20 SetVdsVgs() [2/2] `void SetVdsVgs (`
 `uint32_t Headstage,`
 `int32_t Vds,`
 `int32_t Vgs)`

Sets Vds and Vgs

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vds</i>	Vds in μV
<i>Vgs</i>	Vgs in μV

11.33.3.21 SetVdVs() [1/2] `void SetVdVs (`
`int32_t Vd,`
`int32_t Vs)`

Sets Vd and Vs

Parameters

<i>Vd</i>	Vd in μV
<i>Vs</i>	Vs in μV

11.33.3.22 SetVdVs() [2/2] `void SetVdVs (`
`uint32_t Headstage,`
`int32_t Vd,`
`int32_t Vs)`

Sets Vd and Vs

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vd</i>	Vd in μV
<i>Vs</i>	Vs in μV

11.33.3.23 SetVdVsDAC() [1/2] `void SetVdVsDAC (`
`int16_t Vd,`
`int16_t Vs)`

Sets Vd and Vs

Parameters

<i>Vd</i>	Vd in DAC Units
<i>Vs</i>	Vs in DAC Units

11.33.3.24 SetVdVsDAC() [2/2] `void SetVdVsDAC (`
 `uint32_t Headstage,`
 `int16_t Vd,`
 `int16_t Vs)`

Sets Vd and VS

Parameters

<i>Headstage</i>	The headstage to query.
<i>Vd</i>	Vd in DAC Units
<i>Vs</i>	Vs in DAC Units

11.33.3.25 SetVoltageRange() [1/2] `void SetVoltageRange (`
 `int32_t range)`

Sets the voltage range

Parameters

<i>range</i>	The voltage range in mV
--------------	-------------------------

11.33.3.26 SetVoltageRange() [2/2] `void SetVoltageRange (`
 `uint32_t Headstage,`
 `int32_t range)`

Sets the voltage range

Parameters

<i>Headstage</i>	The headstage to query.
<i>range</i>	The voltage range in mV

11.33.3.27 SetVoltageResolution() [1/2] `void SetVoltageResolution (`
 `int32_t resolution)`

Sets the voltage resolution

Parameters

<i>resolution</i>	The voltage resolution in $\mu\text{V}/\text{digit}$
-------------------	--

11.33.3.28 SetVoltageResolution() [2/2] `void SetVoltageResolution (`
`uint32_t Headstage,`
`int32_t resolution)`

Sets the voltage resolution


Parameters

<i>Headstage</i>	The headstage to query.
<i>resolution</i>	The voltage resolution in $\mu\text{V}/\text{digit}$

11.34 CHiClampDeviceNet Class Reference

[CHiClampDeviceNet](#) is the to control the MCS HiClamp device

Inheritance diagram for CHiClampDeviceNet:

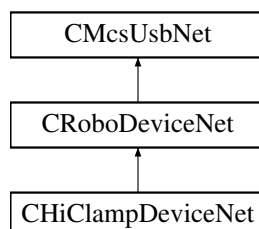

Public Member Functions

- [CHiClampDeviceNet](#) (void)

Properties

- [CRoboDacqNet](#)[^] [RoboDacq](#) [get]

Additional Inherited Members

11.34.1 Detailed Description

[CHiClampDeviceNet](#) is the to control the MCS HiClamp device

11.34.2 Constructor & Destructor Documentation


11.34.2.1 CHiClampDeviceNet() [CHiClampDeviceNet](#) (
`void`)

11.34.3 Property Documentation

11.34.3.1 RoboDacq `CRoboDacqNet`[^] RoboDacq [get]

11.35 CHLADacqNet Class Reference

Inheritance diagram for CHLADacqNet:

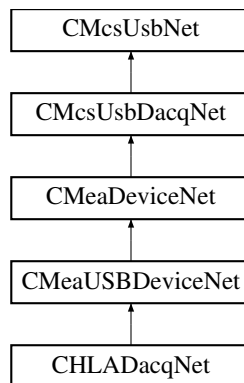

Public Member Functions

- `CHLADacqNet` (void)

Additional Inherited Members


11.35.1 Constructor & Destructor Documentation

11.35.1.1 CHLADacqNet() `CHLADacqNet` (void)

11.36 CHLADeviceNet Class Reference

`CHLADeviceNet` is the to control the MCS HLA device

Inheritance diagram for CHLADeviceNet:

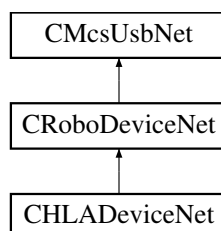

Public Member Functions

- [CHLADeviceNet](#) (void)

Properties

- [CHLADacqNet](#)[^] [HLADacq](#) [get]
- [CSerialPortNet](#)[^] [SerialPort](#) [get]

Additional Inherited Members

11.36.1 Detailed Description

[CHLADeviceNet](#) is the to control the MCS HLA device

11.36.2 Constructor & Destructor Documentation

11.36.2.1 CHLADeviceNet() [CHLADeviceNet](#) (
void)

11.36.3 Property Documentation

11.36.3.1 HLADacq [CHLADacqNet](#)[^] [HLADacq](#) [get]

11.36.3.2 SerialPort [CSerialPortNet](#)[^] [SerialPort](#) [get]

11.37 CMcsUsbDacqNet::CHWInfo Class Reference

Class to provide hardware information about the device.

Classes

- class [CVoltageRangeInfoNet](#)

Public Member Functions

- **CHWInfo** (**CMcsUsbDacqNet**[^] device)
- virtual uint32_t **GetNumberOfHWADCCChannels** ([System::Runtime::InteropServices::Out]int% numberOfChannels)
Get the number of analog channels the device supports.
- virtual uint32_t **GetNumberOfHWDigitalChannels** ([System::Runtime::InteropServices::Out]int% numberOfChannels)
Get the number of digital channels the device supports.
- virtual bool **IsDigitalChannelDedicated** ()
Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.
- virtual uint32_t **GetAvailableSampleRates** ([System::Runtime::InteropServices::Out]System::Collections::Generic::List< int32_t >^% sampleRates)
- virtual System::Collections::Generic::List< int32_t > ^ **GetAvailableVoltageRangesInMicroVolt** (int milliGain)
Gets a List of voltage ranges the device supports.
- virtual System::Collections::Generic::List< **CVoltageRangeInfoNet**[^] > ^ **GetAvailableVoltageRangesInMicroVoltAndStringsInMicroVolt** (int milliGain)
Gets a List of voltage ranges the device supports.

11.37.1 Detailed Description

Class to provide hardware information about the device.

11.37.2 Constructor & Destructor Documentation

11.37.2.1 CHWInfo() **CHWInfo** (
 CMcsUsbDacqNet[^] device)

11.37.3 Member Function Documentation

11.37.3.1 GetAvailableSampleRates() virtual uint32_t GetAvailableSampleRates (
 [System::Runtime::InteropServices::Out] System::Collections::Generic::List< int32_t >^% sampleRates) [virtual]

11.37.3.2 GetAvailableVoltageRangesInMicroVolt() virtual System::Collections::Generic::List<int32_t> ^ GetAvailableVoltageRangesInMicroVolt (
 int milliGain) [virtual]

Gets a List of voltage ranges the device supports.

The List is scaled by the gain factor given as argument to this function. Use "1000" as scale factor for backwards compatibility. To get a list of voltage ranges for the headstage, obtain the gain of the headstage with the **GetGain()** call and use the result in the milliGain parameter. To get a list of voltage ranges for the analog inputs of the interfaceboard, obtain the gain of the analog inputs with the **GetAnalogGain()** call and use the result in the milliGain parameter.

Parameters

<i>milliGain</i>	The gain factor (in milliGain) used to scale the list of voltage ranges.
------------------	--

Returns

List of voltage ranges in μV .

11.37.3.3 GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt() `virtual System::Collections::Generic::List<CVoltageRangeInfoNet^> ^ GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt (int milliGain) [virtual]`

Gets a List of voltage ranges the device supports.

The List is scaled by the gain factor given as argument to this function. Use "1000" as scale factor for backwards compatibility. Each list entry contains the voltage range as an integer and as a string. To get a list of voltage ranges for the headstage, obtain the gain of the headstage with the GetGain() call and use the result in the milliGain parameter. To get a list of voltage ranges for the analog inputs of the interfaceboard, obtain the gain of the analog inputs with the GetAnalogGain() call and use the result in the milliGain parameter.

Parameters

<i>milliGain</i>	The gain factor (in milliGain) used to scale the list of voltage ranges.
------------------	--

Returns

List of voltage ranges in μV .

11.37.3.4 GetNumberOfHWADCCChannels() `virtual uint32_t GetNumberOfHWADCCChannels ([System::Runtime::InteropServices::Out] int% numberOfChannels) [virtual]`

Get the number of analog channels the device supports.

Parameters

<i>numberOfChannels</i>	Number of analog channels the device supports.
-------------------------	--

Returns

Error Status. 0 on success.

11.37.3.5 GetNumberOfHWDigitalChannels() `virtual uint32_t GetNumberOfHWDigitalChannels ([System::Runtime::InteropServices::Out] int% numberOfChannels) [virtual]`

Get the number of digital channels the device supports.

Parameters

<i>numberOfChannels</i>	Number of digital channels the device supports.
-------------------------	---

Returns

Error Status. 0 on success.

11.37.3.6 IsDigitalChannelDedicated() `virtual bool IsDigitalChannelDedicated () [virtual]`


Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

Returns

false when the digital channel replaces an analog channel when enabled, true when the digital channels is appended to the analog channels when enabled.

11.38 CIntanMea_FunctionNet Class Reference

Inheritance diagram for CIntanMea_FunctionNet:

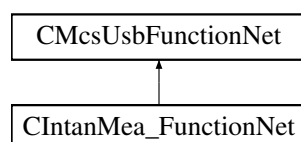

Public Member Functions

- [CIntanMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] intalMea_FunctionPointerContainer)
- [CIntanMea_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- [int](#) [GetUpperFrequencyByIndex](#) (unsigned short index)
- [int](#) [GetLowerFrequencyByIndex](#) (unsigned short index)
- [int64_t](#) [GetDSPHighPassByIndex](#) (unsigned short index)
- [int](#) [GetIntanRegister](#) (unsigned short chip, unsigned short registernumber)
- [int](#) [GetImpedanceResult](#) (unsigned short channel)
- [void](#) [SetBandwidthByIndex](#) (int upper_index, int lower_index)
- [void](#) [SetDSPHighPassByIndex](#) (int index)
- [void](#) [AmplifierSettle](#) ()
- [void](#) [SetIntanRegister](#) (unsigned short register_number, int value)
- [void](#) [SetDiagnosticMode](#) (unsigned char onoff)
- [void](#) [BeginImpedanceCheck](#) (array< int >[^] config_values)

Additional Inherited Members

11.38.1 Constructor & Destructor Documentation

11.38.1.1 CIntanMea_FunctionNet() [1/2] `CIntanMea_FunctionNet (`
`CMcsUsbNet^ mcsusb,`
`CMcsUsbFunctionPointerContainer^ intalMea_FunctionPointerContainer)`

11.38.1.2 CIntanMea_FunctionNet() [2/2] `CIntanMea_FunctionNet (`
`CMcsUsbNet^ mcsusb)`

11.38.2 Member Function Documentation

11.38.2.1 AmplifierSettle() `void AmplifierSettle ()`

11.38.2.2 BeginImpedanceCheck() `void BeginImpedanceCheck (`
`array< int >^ config_values)`

11.38.2.3 GetDSPHighPassByIndex() `int64_t GetDSPHighPassByIndex (`
`unsigned short index)`

11.38.2.4 GetImpedanceResult() `int GetImpedanceResult (`
`unsigned short channel)`

11.38.2.5 GetIntanRegister() `int GetIntanRegister (`
`unsigned short chip,`
`unsigned short registernumber)`

11.38.2.6 GetLowerFrequencyByIndex() `int GetLowerFrequencyByIndex (`
`unsigned short index)`

11.38.2.7 GetUpperFrequencyByIndex() `int GetUpperFrequencyByIndex (`
`unsigned short index)`

11.38.2.8 SetBandwidthByIndex() `void SetBandwidthByIndex (`
`int upper_index,`
`int lower_index)`

11.38.2.9 SetDiagnosticMode() `void SetDiagnosticMode (`
`unsigned char onoff)`


11.38.2.10 SetDSPHighPassByIndex() `void SetDSPHighPassByIndex (`
`int index)`

11.38.2.11 SetIntanRegister() `void SetIntanRegister (`
`unsigned short register_number,`
`int value)`

11.39 CInterfaceboard2FunctionNet Class Reference

[CInterfaceboard2FunctionNet](#) is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboard2FunctionNet:

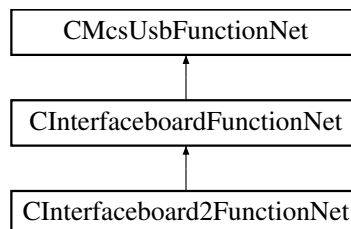

Public Member Functions

- [CInterfaceboard2FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pInterfaceboard2FunctionPointerContainer)
Initializes a new instance of the [CInterfaceboard2FunctionNet](#) class.
- [CInterfaceboard2FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CInterfaceboard2FunctionNet](#) ()
- [!CInterfaceboard2FunctionNet](#) ()
- void [SetIoVoltage](#) ([IoVoltageEnumNet](#) ioVoltage)
Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.
- [IoVoltageEnumNet](#) [GetIoVoltage](#) ()
Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

Additional Inherited Members

11.39.1 Detailed Description

[CInterfaceboard2FunctionNet](#) is the class to control the Interfaceboard

11.39.2 Constructor & Destructor Documentation

11.39.2.1 CInterfaceboard2FunctionNet() [1/2] [CInterfaceboard2FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pInterfaceboard2FunctionPointerContainer)

Initializes a new instance of the [CInterfaceboard2FunctionNet](#) class.

11.39.2.2 CInterfaceboard2FunctionNet() [2/2] [CInterfaceboard2FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

11.39.2.3 ~CInterfaceboard2FunctionNet() virtual [~CInterfaceboard2FunctionNet](#) () [virtual]

11.39.2.4 !CInterfaceboard2FunctionNet() [!CInterfaceboard2FunctionNet](#) ()

11.39.3 Member Function Documentation

11.39.3.1 GetIoVoltage() `IoVoltageEnumNet GetIoVoltage ()`

Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

Returns

Enum for the IO Voltage (3.3V or 5.0V).

11.39.3.2 SetIoVoltage() `void SetIoVoltage (IoVoltageEnumNet ioVoltage)`

Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.


Parameters

<code>ioVoltage</code>	Enum for the I/O Voltage (3.3V or 5.0V).
------------------------	--

11.40 CInterfaceboardFunctionNet Class Reference

[CInterfaceboardFunctionNet](#) is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboardFunctionNet:

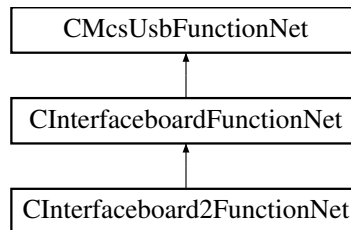

Public Member Functions

- [CInterfaceboardFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pInterfaceboard↔
FunctionPointerContainer)
Initializes a new instance of the [CInterfaceboardFunctionNet](#) class.
- [CInterfaceboardFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CInterfaceboardFunctionNet](#) ()
- [!CInterfaceboardFunctionNet](#) ()
- void [SetCardinalDacqSamplerate](#) (uint32_t samplerate)
Sets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz
- uint32_t [GetCardinalDacqSamplerate](#) ()
Gets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz
- void [SetCardinalStgOutputrate](#) (uint32_t outputrate)
Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz
- uint32_t [GetCardinalStgOutputrate](#) ()
Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Additional Inherited Members

11.40.1 Detailed Description

[CInterfaceboardFunctionNet](#) is the class to control the Interfaceboard

11.40.2 Constructor & Destructor Documentation

11.40.2.1 CInterfaceboardFunctionNet() [1/2] [CInterfaceboardFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pInterfaceboardFunctionPointerContainer*)

Initializes a new instance of the [CInterfaceboardFunctionNet](#) class.

11.40.2.2 CInterfaceboardFunctionNet() [2/2] [CInterfaceboardFunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

11.40.2.3 ~CInterfaceboardFunctionNet() virtual [~CInterfaceboardFunctionNet](#) () [virtual]

11.40.2.4 "!CInterfaceboardFunctionNet() [!CInterfaceboardFunctionNet](#) ()

11.40.3 Member Function Documentation

11.40.3.1 GetCardinalDacqSamplerate() uint32_t GetCardinalDacqSamplerate ()

Gets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz

Returns

The samplerate in Hz.

11.40.3.2 GetCardinalStgOutputrate() uint32_t GetCardinalStgOutputrate ()

Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Returns

The output rate in Hz.

11.40.3.3 SetCardinalDacqSamplerate() void SetCardinalDacqSamplerate (
 uint32_t *samplerate*)

Sets the fundamental/cardinal data acquisition samplerate of the Interfaceboard, default is 50 kHz

Parameters

<i>samplerate</i>	The samplerate in Hz.
-------------------	-----------------------

11.40.3.4 SetCardinalStgOutputrate() `void SetCardinalStgOutputrate (`
`uint32_t outputrate)`

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz


Parameters

<i>outputrate</i>	The output rate in Hz.
-------------------	------------------------

11.41 CLIH3DeviceNet Class Reference

[CLIH3DeviceNet](#) is the class to access the HEKA LIH3 device.

Inheritance diagram for CLIH3DeviceNet:

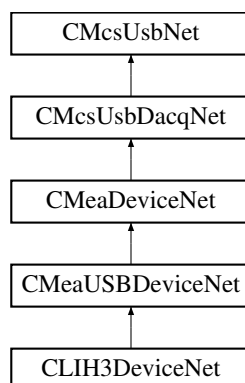

Public Member Functions

- [CLIH3DeviceNet](#) ()
Initializes a new instance of the [CLIH3DeviceNet](#) class.
- virtual [~CLIH3DeviceNet](#) ()
- [!CLIH3DeviceNet](#) ()
- void [DummyCommand](#) (uint32_t dummyParameter)
Dummy command to show how to use the DLL.
- void [SetEEpromPage](#) (uint32_t EEpromStartAddress, array< int8_t >^ EEpromData, [LIH30_EPC10_Bus_EnumNet](#) epc10bus)
Writes into EEprom on the EPC10 EEPROM
- array< int8_t >^ [GetEEpromPage](#) (uint32_t EEpromStartAddress, int EEpromData_Length, [LIH30_EPC10_Bus_EnumNet](#) epc10bus)
Reads the requested amount of EEprom byte from the EPC10 EEPROM
- void [SetSampleInterval](#) (uint32_t SampleInterval)

- Sets the Sample Interval for the DACQ and Stimulation*

 - uint32_t [GetSampleInterval](#) ()
- Gets the Sample Interval for the DACQ and Stimulation*

 - void [SetAdcOffset](#) (LIH30_ADC_Channel_EnumNet AdcChannel, int32_t Offset)
- Sets the ADC offset of the DACQ for a single channel*

 - int32_t [GetAdcOffset](#) (LIH30_ADC_Channel_EnumNet AdcChannel)
- Gets the ADC offset of the DACQ for a single channel*

 - void [SetAdcOffsetPermanent](#) (LIH30_ADC_Channel_EnumNet AdcChannel)
- Writes the ADC offset of the DACQ for a single channel to permanent EEPROM memory*

 - void [ErasePermanentAdcOffset](#) (LIH30_ADC_Channel_EnumNet AdcChannel)
- Deletes the ADC offset of the DACQ for a single channel in permanent EEPROM memory*

 - uint32_t [ReadClipping](#) (LIH30_EPC10_Bus_EnumNet epc10bus)
- Gets the clipping information*

 - void [SetDigOutState](#) (uint16_t DigOutState)
- Writes to the LIH30 digital output*

 - uint16_t [GetDigInState](#) ()
- Reads from the LIH30 digital input*

 - void [SendCommand](#) (LIH30_EPC10_Bus_EnumNet epc10bus, uint16_t Command)
- Send command to the EPC10*

 - uint16_t [GetDacqRunStatus](#) ()
- Gets the data acquisition running status*

 - void [SetDacUseldleValue](#) (uint32_t DacChannel, bool Useldle)
- Sets if the DAC Idle value is used after stimulation*

 - bool [GetDacUseldleValue](#) (uint32_t DacChannel)
- Gets if the DAC Idle value is used after stimulation*

 - void [SetDacIdleValue](#) (uint32_t DacChannel, int32_t IdleValue)
- Sets the DAC Idle value*

 - int32_t [GetDacIdleValue](#) (uint32_t DacChannel)
- Gets the DAC Idle value*

 - void [EnableUserTrigger](#) (bool enable)
- Enables the User Trigger*

 - bool [IsUserTriggerEnabled](#) ()
- Is the User Trigger enabled*

 - void [SetDacOffset](#) (LIH30_DAC_Channel_EnumNet DacChannel, int32_t Offset)
- Sets the offset of a DAC channel.*

 - int32_t [GetDacOffset](#) (LIH30_DAC_Channel_EnumNet DacChannel)
- Gets the offset of a DAC channel.*

 - void [SetDacOffsetPermanent](#) (LIH30_DAC_Channel_EnumNet DacChannel)
- Writes the DAC offset of the STG for a single channel to permanent EEPROM memory*

 - void [ErasePermanentDacOffset](#) (LIH30_DAC_Channel_EnumNet DacChannel)
- Deletes the DAC offset of the STG for a single channel in permanent EEPROM memory*

 - void [SetAudioOutDacParameter](#) (uint32_t Frequency, uint32_t Amplification)
- Sets the parameter of the audio DAC output.*

 - void [GetAudioOutDacParameter](#) ([System::Runtime::InteropServices::Out]uint32_t% Frequency, [System::Runtime::InteropServices::Out]uint32_t% Amplification)
- Gets the parameter of the audio DAC output.*

 - String ^ [ReadUARTData](#) ()
- Reads the config string from the device.*

 - void [WriteUARTData](#) (String^ commandString)
- Write the command string to the device.*

Properties

- [CStimulusFunctionNet](#)[^] [StimulusFunction](#) [get]

Additional Inherited Members

11.41.1 Detailed Description

[CLIH3DeviceNet](#) is the class to access the HEKA LIH3 device.

11.41.2 Constructor & Destructor Documentation

11.41.2.1 [CLIH3DeviceNet\(\)](#) [CLIH3DeviceNet](#) ()

Initializes a new instance of the [CLIH3DeviceNet](#) class.

11.41.2.2 [~CLIH3DeviceNet\(\)](#) [virtual ~CLIH3DeviceNet](#) () [virtual]

11.41.2.3 ["!CLIH3DeviceNet\(\)](#) [!CLIH3DeviceNet](#) ()

11.41.3 Member Function Documentation

11.41.3.1 [DummyCommand\(\)](#) [void DummyCommand](#) ([uint32_t dummyParameter](#))

Dummy command to show how to use the DLL.

Parameters

<i>dummyParameter</i>	parameter to send to the device
-----------------------	---------------------------------

11.41.3.2 [EnableUserTrigger\(\)](#) [void EnableUserTrigger](#) ([bool enable](#))

Enables the User Trigger

Parameters

<i>enable</i>	Enable
---------------	--------

11.41.3.3 ErasePermanentAdcOffset() `void ErasePermanentAdcOffset (
 LIH30_ADC_Channel_EnumNet AdcChannel)`

Deletes the ADC offset of the DACQ for a single channel in permanent EEPROM memory

Parameters

<i>AdcChannel</i>	The ADC channel
-------------------	-----------------

11.41.3.4 ErasePermanentDacOffset() `void ErasePermanentDacOffset (
 LIH30_DAC_Channel_EnumNet DacChannel)`

Deletes the DAC offset of the STG for a single channel in permanent EEPROM memory

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

11.41.3.5 GetAdcOffset() `int32_t GetAdcOffset (
 LIH30_ADC_Channel_EnumNet AdcChannel)`

Gets the ADC offset of the DACQ for a single channel

Parameters

<i>AdcChannel</i>	The ADC channel
-------------------	-----------------

Returns

The offset for the given channel number

11.41.3.6 GetAudioOutDacParameter() `void GetAudioOutDacParameter (
 [System::Runtime::InteropServices::Out] uint32_t% Frequency,
 [System::Runtime::InteropServices::Out] uint32_t% Amplification)`

Gets the parameter of the audio DAC output.

Parameters

<i>Frequency</i>	Frequency(1 - 25000 Hz)
<i>Amplification</i>	Amplification(0 - 0xFFFF)

11.41.3.7 GetDacIdleValue() `int32_t GetDacIdleValue (`
`uint32_t DacChannel)`

Gets the DAC Idle value

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

Returns

The idle value

11.41.3.8 GetDacOffset() `int32_t GetDacOffset (`
`LIH30_DAC_Channel_EnumNet DacChannel)`

Gets the offset of a DAC channel.

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

Returns

The offset for the given channel number

11.41.3.9 GetDacqRunStatus() `uint16_t GetDacqRunStatus ()`

Gets the data acquisition running status

Returns

The status (1: running / 0: stopped)

11.41.3.10 GetDacUseIdleValue() `bool GetDacUseIdleValue (`
`uint32_t DacChannel)`

Gets if the DAC Idle value is used after stimulation

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

Returns

Use idle value

11.41.3.11 GetDigInState() `uint16_t GetDigInState ()`

Reads from the LIH30 digital input

Returns

The bit mask defining the digital input state

11.41.3.12 GetEEpromPage() `array<int8_t> ^ GetEEpromPage (uint32_t EEpromStartAddress, int EEpromData_Length, LIH30_EPC10_Bus_EnumNet epc10bus)`

Reads the requested amount of EEprom byte from the EPC10 EEPROM

Parameters

<i>EEpromStartAddress</i>	start address of memory area to read from
<i>EEpromData_Length</i>	The maximal length of EEpromData.
<i>epc10bus</i>	The EPC10 bus

Returns

pointer to internal memory for the requested amount of data

11.41.3.13 GetSampleInterval() `uint32_t GetSampleInterval ()`

Gets the Sample Interval for the DACQ and Stimulation

Returns

Sample Interval configured on the device

11.41.3.14 IsUserTriggerEnabled() `bool IsUserTriggerEnabled ()`

Is the User Trigger enabled

Returns

Enabled

11.41.3.15 ReadClipping() `uint32_t ReadClipping (`
`LIH30_EPC10_Bus_EnumNet epc10bus)`

Gets the clipping information

Parameters

<i>epc10bus</i>	The EPC10 bus
-----------------	---------------

Returns

The clipping value

11.41.3.16 ReadUARTData() `String ^ ReadUARTData ()`

Reads the config string from the device.

Returns

The config string.

11.41.3.17 SendCommand() `void SendCommand (`
`LIH30_EPC10_Bus_EnumNet epc10bus,`
`uint16_t Command)`

Send command to the EPC10

Parameters

<i>epc10bus</i>	The EPC10 bus
<i>Command</i>	The command

11.41.3.18 SetAdcOffset() `void SetAdcOffset (`
 `LIH30_ADC_Channel_EnumNet AdcChannel,`
 `int32_t Offset)`

Sets the ADC offset of the DACQ for a single channel

Parameters

<i>AdcChannel</i>	The ADC channel
<i>Offset</i>	The offset for the given channel number

11.41.3.19 SetAdcOffsetPermanent() `void SetAdcOffsetPermanent (`
 `LIH30_ADC_Channel_EnumNet AdcChannel)`

Writes the ADC offset of the DACQ for a single channel to permanent EEPROM memory

Parameters

<i>AdcChannel</i>	The ADC channel
-------------------	-----------------

11.41.3.20 SetAudioOutDacParameter() `void SetAudioOutDacParameter (`
 `uint32_t Frequency,`
 `uint32_t Amplification)`

Sets the parameter of the audio DAC output.

Parameters

<i>Frequency</i>	Frequency(1 - 25000 Hz)
<i>Amplification</i>	Amplification(0 - 0xFFFF)

11.41.3.21 SetDacIdleValue() `void SetDacIdleValue (`
 `uint32_t DacChannel,`
 `int32_t IdleValue)`

Sets the DAC Idle value

Parameters

<i>DacChannel</i>	The DAC channel
<i>IdleValue</i>	The idle value

11.41.3.22 SetDacOffset() `void SetDacOffset (`
 `LIH30_DAC_Channel_EnumNet DacChannel,`
 `int32_t Offset)`

Sets the offset of a DAC channel.

Parameters

<i>DacChannel</i>	The DAC channel
<i>Offset</i>	The offset for the given channel number

11.41.3.23 SetDacOffsetPermanent() `void SetDacOffsetPermanent (`
 `LIH30_DAC_Channel_EnumNet DacChannel)`

Writes the DAC offset of the STG for a single channel to permanent EEPROM memory

Parameters

<i>DacChannel</i>	The DAC channel
-------------------	-----------------

11.41.3.24 SetDacUseIdleValue() `void SetDacUseIdleValue (`
 `uint32_t DacChannel,`
 `bool UseIdle)`

Sets if the DAC Idle value is used after stimulation

Parameters

<i>DacChannel</i>	The DAC channel
<i>UseIdle</i>	Use idle value

11.41.3.25 SetDigOutState() `void SetDigOutState (`
 `uint16_t DigOutState)`

Writes to the LIH30 digital output

Parameters

<i>DigOutState</i>	The bit mask defining the digital output state
--------------------	--

11.41.3.26 SetEepromPage() `void SetEepromPage (`

```
uint32_t EepromStartAddress,
array< int8_t >^ EepromData,
LIH30_EPC10_Bus_EnumNet epc10bus )
```

Writes into EEprom on the EPC10 EEPROM

Parameters

<i>EepromStartAddress</i>	start address of memory area to write to
<i>EepromData</i>	pointer to internal memory for the supported amount of data
<i>epc10bus</i>	The EPC10 bus

11.41.3.27 SetSampleInterval() void SetSampleInterval (
uint32_t *SampleInterval*)

Sets the Sample Interval for the DACQ and Stimulation

Parameters

<i>SampleInterval</i>	between the samples, Sample interval is available from 1 to 4194303
-----------------------	---

11.41.3.28 WriteUARTData() void WriteUARTData (
String^ *commandString*)

Write the command string to the device.

Parameters


<i>commandString</i>	The config string.
----------------------	--------------------

11.41.4 Property Documentation

11.41.4.1 StimulusFunction CStimulusFunctionNet^ StimulusFunction [get]

11.42 CMcsBus_AxisParametersNet Class Reference

Inheritance diagram for CMcsBus_AxisParametersNet:

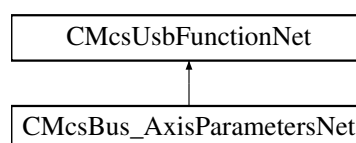

Public Member Functions

- [CMcsBus_AxisParametersNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_AxisParametersNet](#) (void)
- void [SetAxisParametersEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, unsigned int parameter)
- void [SetAxisParametersEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, int parameter)
- unsigned int [GetAxisParametersUnsignedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)
- int [GetAxisParametersSignedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)

Additional Inherited Members

11.42.1 Constructor & Destructor Documentation

11.42.1.1 [CMcsBus_AxisParametersNet\(\)](#) [CMcsBus_AxisParametersNet](#) (
 [CMcsUsbNet](#)[^] device)

11.42.1.2 [~CMcsBus_AxisParametersNet\(\)](#) [~CMcsBus_AxisParametersNet](#) (
 void)

11.42.2 Member Function Documentation

11.42.2.1 [GetAxisParametersSignedEeprom\(\)](#) int [GetAxisParametersSignedEeprom](#) (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *index*)


11.42.2.2 [GetAxisParametersUnsignedEeprom\(\)](#) unsigned int [GetAxisParametersUnsignedEeprom](#) (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *index*)

11.42.2.3 SetAxisParametersEeprom() [1/2] `void SetAxisParametersEeprom (`
`unsigned char busnumber,`
`unsigned char busaddress,`
`unsigned char axis,`
`unsigned short index,`
`int parameter)`

11.42.2.4 SetAxisParametersEeprom() [2/2] `void SetAxisParametersEeprom (`
`unsigned char busnumber,`
`unsigned char busaddress,`
`unsigned char axis,`
`unsigned short index,`
`unsigned int parameter)`

11.43 CMcsBus_ExtensionNet Class Reference

Inheritance diagram for CMcsBus_ExtensionNet:

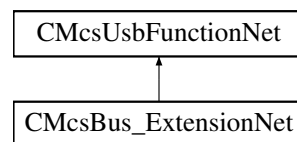

Public Member Functions

- [CMcsBus_ExtensionNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_ExtensionNet](#) (void)
- void [SetLEDSwitch](#) (unsigned char busnumber, unsigned char busaddress, unsigned short LEDSwitch)
- unsigned short [GetLEDSwitch](#) (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.43.1 Constructor & Destructor Documentation

11.43.1.1 CMcsBus_ExtensionNet() [CMcsBus_ExtensionNet](#) (
[CMcsUsbNet](#)[^] device)

11.43.1.2 ~CMcsBus_ExtensionNet() [~CMcsBus_ExtensionNet](#) (
void)


11.43.2 Member Function Documentation

11.43.2.1 GetLEDSwitch() unsigned short GetLEDSwitch (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.43.2.2 SetLEDSwitch() void SetLEDSwitch (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *LEDSwitch*)

11.44 CMcsBus_FYIExtensionNet Class Reference

Inheritance diagram for CMcsBus_FYIExtensionNet:

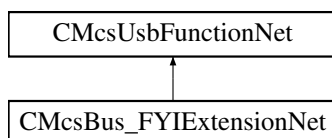

Public Member Functions

- `CMcsBus_FYIExtensionNet` (`CMcsUsbNet`[^] device)
- `~CMcsBus_FYIExtensionNet` (void)
- void `SetValves` (unsigned char busnumber, unsigned char busaddress, unsigned int states)
- unsigned int `GetValves` (unsigned char busnumber, unsigned char busaddress)
- void `SetDIO` (unsigned char busnumber, unsigned char busaddress, unsigned short io)
- unsigned short `GetDIO` (unsigned char busnumber, unsigned char busaddress)
- void `SetSingleHeater` (unsigned char busnumber, unsigned char busaddress, short index, unsigned short power)
- unsigned short `GetSingleHeater` (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members

11.44.1 Constructor & Destructor Documentation

11.44.1.1 CMcsBus_FYIExtensionNet() `CMcsBus_FYIExtensionNet` (
 `CMcsUsbNet`[^] *device*)

11.44.1.2 `~CMcsBus_FYIExtensionNet()` `~CMcsBus_FYIExtensionNet` (
 `void`)

11.44.2 Member Function Documentation

11.44.2.1 `GetDIO()` `unsigned short GetDIO` (
 `unsigned char busnumber`,
 `unsigned char busaddress`)

11.44.2.2 `GetSingleHeater()` `unsigned short GetSingleHeater` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `short index`)

11.44.2.3 `GetValves()` `unsigned int GetValves` (
 `unsigned char busnumber`,
 `unsigned char busaddress`)


11.44.2.4 `SetDIO()` `void SetDIO` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `unsigned short io`)

11.44.2.5 `SetSingleHeater()` `void SetSingleHeater` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `short index`,
 `unsigned short power`)

11.44.2.6 `SetValves()` `void SetValves` (
 `unsigned char busnumber`,
 `unsigned char busaddress`,
 `unsigned int states`)

11.45 CMcsBus_MotorControlNet Class Reference

Inheritance diagram for CMcsBus_MotorControlNet:

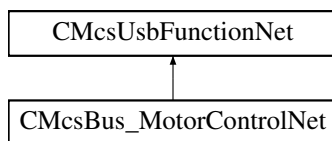

Public Member Functions

- [CMcsBus_MotorControlNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_MotorControlNet](#) (void)
- void [SetMCScalingFactorEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int [GetMCScalingFactorEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCScalingFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int [GetMCScalingFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCMaxSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCMaxSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravelEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravelEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCMaxCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCMaxCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRegulatorGainEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short [GetMCRegulatorGainEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRegulatorGain](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short [GetMCRegulatorGain](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCMaxAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void [SetMCMaxAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCMaxAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short [GetMCStandbyCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short [GetMCStandbyCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyTimeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short [GetMCStandbyTimeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCStandbyTime](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short [GetMCStandbyTime](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCBreakCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCBreakCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCBreakCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCBreakCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCConfigEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short [GetMCConfigEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCConfig](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short [GetMCConfig](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short [GetMCSpeedEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short [GetMCSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAccelerationEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAcceleration](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReferenceCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCReferenceCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReferenceCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCReferenceCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [RoboCurrentModeEnumNet](#) mode)
- [RoboCurrentModeEnumNet](#) [GetMCCurrentModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void [SetMCCurrentMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [RoboCurrentModeEnumNet](#) mode)
- [RoboCurrentModeEnumNet GetMCCurrentMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAxisRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short revision)
- unsigned short [GetMCAxisRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedUnitEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int32_t speedunit)
- int32_t [GetMCSpeedUnitEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, bool OnOff_status)
- bool [GetMCOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCSpeedShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short [GetMCSpeedShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCAccelerationShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short [GetMCAccelerationShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short [GetMCCurrentShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCMaxTravelShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int [GetMCMaxTravelShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int [GetMCCurrentPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCNewPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int [GetMCNewPosition](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- short [GetMCCurrentSpeed](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [StartMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCRotation](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char onoff)
- unsigned short [GetMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCReference](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char switch_enable, unsigned char switch_polarity)
- unsigned char [GetMCReference](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned char% switch_port)
- void [StopMCMovement](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetMCCurrentModeShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [RoboCurrentModeEnumNet](#) mode)
- [RoboCurrentModeEnumNet GetMCCurrentModeShortCommand](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short [GetMCPhase](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short [GetMCPhaseOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void [SetSubChannel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short subchannel)
- unsigned short [GetSubChannel](#) (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

Additional Inherited Members

11.45.1 Constructor & Destructor Documentation

11.45.1.1 CMcsBus_MotorControlNet() `CMcsBus_MotorControlNet (CMcsUsbNet^ device)`

11.45.1.2 ~CMcsBus_MotorControlNet() `~CMcsBus_MotorControlNet (void)`

11.45.2 Member Function Documentation

11.45.2.1 GetMCAcceleration() `unsigned short GetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

11.45.2.2 GetMCAccelerationEeprom() `unsigned short GetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

11.45.2.3 GetMCAccelerationShortCommand() `unsigned short GetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

11.45.2.4 GetMCAXisRevisionEeprom() `unsigned short GetMCAXisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)`

11.45.2.5 GetMCBreakCurrent() short GetMCBreakCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.6 GetMCBreakCurrentEeprom() short GetMCBreakCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.7 GetMCConfig() unsigned short GetMCConfig (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.8 GetMCConfigEeprom() unsigned short GetMCConfigEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.9 GetMCCurrent() short GetMCCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.10 GetMCCurrentEeprom() short GetMCCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.11 GetMCCurrentMode() [RoboCurrentModeEnumNet](#) GetMCCurrentMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.12 GetMCCurrentModeEeprom() [RoboCurrentModeEnumNet](#) GetMCCurrentModeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.13 GetMCCurrentModeShortCommand() [RoboCurrentModeEnumNet](#) GetMCCurrentModeShort↵
Command (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.14 GetMCCurrentPosition() int GetMCCurrentPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.15 GetMCCurrentShortCommand() short GetMCCurrentShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.16 GetMCCurrentSpeed() short GetMCCurrentSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.17 GetMCMaxAcceleration() unsigned short GetMCMaxAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.18 GetMCMaxAccelerationEeprom() unsigned short GetMCMaxAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.19 GetMCMaxCurrent() short GetMCMaxCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.20 GetMCMaxCurrentEeprom() short GetMCMaxCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.21 GetMCMaxSpeed() unsigned short GetMCMaxSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.22 GetMCMaxSpeedEeprom() unsigned short GetMCMaxSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.23 GetMCMaxTravel() int GetMCMaxTravel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.24 GetMCMaxTravelEeprom() int GetMCMaxTravelEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.25 GetMCMaxTravelShortCommand() int GetMCMaxTravelShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

- 11.45.2.26 GetMCMovement()** unsigned short GetMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 11.45.2.27 GetMCNewPosition()** int GetMCNewPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 11.45.2.28 GetMCOutputOnOff()** bool GetMCOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 11.45.2.29 GetMCPhase()** unsigned short GetMCPhase (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 11.45.2.30 GetMCPhaseOffset()** unsigned short GetMCPhaseOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)
- 11.45.2.31 GetMCReference()** unsigned char GetMCReference (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned char% *switch_port*)
- 11.45.2.32 GetMCReferenceCurrent()** short GetMCReferenceCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.33 GetMCReferenceCurrentEeprom() short GetMCReferenceCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.34 GetMCRegulatorGain() short GetMCRegulatorGain (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.35 GetMCRegulatorGainEeprom() short GetMCRegulatorGainEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.36 GetMCScalingFactor() int GetMCScalingFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.37 GetMCScalingFactorEeprom() int GetMCScalingFactorEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.38 GetMCSpeed() short GetMCSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.39 GetMCSpeedEeprom() unsigned short GetMCSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.40 GetMCSpeedShortCommand() short GetMCSpeedShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.41 GetMCSpeedUnitEeprom() int32_t GetMCSpeedUnitEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.42 GetMCStandbyCurrent() short GetMCStandbyCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.43 GetMCStandbyCurrentEeprom() short GetMCStandbyCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.44 GetMCStandbyTime() short GetMCStandbyTime (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.45 GetMCStandbyTimeEeprom() short GetMCStandbyTimeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.46 GetSubChannel() unsigned short GetSubChannel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.47 SetMCAcceleration() void SetMCAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

11.45.2.48 SetMCAccelerationEeprom() void SetMCAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

11.45.2.49 SetMCAccelerationShortCommand() void SetMCAccelerationShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

11.45.2.50 SetMCAxisRevisionEeprom() void SetMCAxisRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *revision*)

11.45.2.51 SetMCBreakCurrent() void SetMCBreakCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.52 SetMCBreakCurrentEeprom() void SetMCBreakCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.53 SetMCConfig() void SetMCConfig (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *config*)

11.45.2.54 SetMCConfigEeprom() void SetMCConfigEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *config*)

11.45.2.55 SetMCCurrent() void SetMCCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.56 SetMCCurrentEeprom() void SetMCCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.57 SetMCCurrentMode() void SetMCCurrentMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

11.45.2.58 SetMCCurrentModeEeprom() void SetMCCurrentModeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

11.45.2.59 SetMCCurrentModeShortCommand() void SetMCCurrentModeShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 RoboCurrentModeEnumNet *mode*)

11.45.2.60 SetMCCurrentPosition() void SetMCCurrentPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *position*)

11.45.2.61 SetMCCurrentShortCommand() void SetMCCurrentShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.62 SetMCMaxAcceleration() void SetMCMaxAcceleration (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

11.45.2.63 SetMCMaxAccelerationEeprom() void SetMCMaxAccelerationEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *acceleration*)

11.45.2.64 SetMCMaxCurrent() void SetMCMaxCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.65 SetMCMaxCurrentEeprom() void SetMCMaxCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.66 SetMCMaxSpeed() void SetMCMaxSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

11.45.2.67 SetMCMaxSpeedEeprom() void SetMCMaxSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

11.45.2.68 SetMCMaxTravel() void SetMCMaxTravel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

11.45.2.69 SetMCMaxTravelEeprom() void SetMCMaxTravelEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

11.45.2.70 SetMCMaxTravelShortCommand() void SetMCMaxTravelShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *travel*)

11.45.2.71 SetMCNewPosition() void SetMCNewPosition (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *position*)

11.45.2.72 SetMCOutputOnOff() void SetMCOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 bool *OnOff_status*)

11.45.2.73 SetMCReference() void SetMCReference (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned char *switch_enable*,
 unsigned char *switch_polarity*)

11.45.2.74 SetMCReferenceCurrent() void SetMCReferenceCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.75 SetMCReferenceCurrentEeprom() void SetMCReferenceCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *current*)

11.45.2.76 SetMCRegulatorGain() void SetMCRegulatorGain (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *gain*)

11.45.2.77 SetMCRegulatorGainEeprom() void SetMCRegulatorGainEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *gain*)

11.45.2.78 SetMCRotation() void SetMCRotation (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned char *onoff*)

11.45.2.79 SetMCScalingFactor() void SetMCScalingFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *factor*)

11.45.2.80 SetMCScalingFactorEeprom() void SetMCScalingFactorEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *factor*)

11.45.2.81 SetMCSpeed() void SetMCSpeed (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *speed*)

11.45.2.82 SetMCSpeedEeprom() void SetMCSpeedEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *speed*)

11.45.2.83 SetMCSpeedShortCommand() void SetMCSpeedShortCommand (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *speed*)

11.45.2.84 SetMCSpeedUnitEeprom() void SetMCSpeedUnitEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int32_t *speedunit*)

11.45.2.85 SetMCStandbyCurrent() void SetMCStandbyCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *percent*)

11.45.2.86 SetMCStandbyCurrentEeprom() void SetMCStandbyCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *percent*)

11.45.2.87 SetMCStandbyTime() void SetMCStandbyTime (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *t*)

11.45.2.88 SetMCStandbyTimeEeprom() void SetMCStandbyTimeEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *t*)


11.45.2.89 SetSubChannel() void SetSubChannel (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *subchannel*)

11.45.2.90 StartMCMovement() void StartMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.45.2.91 StopMCMovement() void StopMCMovement (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*)

11.46 CMcsBus_SensorNet Class Reference

Inheritance diagram for CMcsBus_SensorNet:

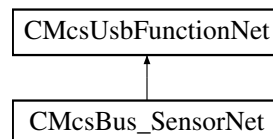

Public Member Functions

- [CMcsBus_SensorNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_SensorNet](#) (void)
- void [SetMinimalThreshold](#) (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short [GetMinimalThreshold](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetDetectionThreshold](#) (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short [GetDetectionThreshold](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetLatency](#) (unsigned char busnumber, unsigned char busaddress, unsigned short latency)
- unsigned short [GetLatency](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetBubbleStatus](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetLatencyCounter](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetDetectorValue](#) (unsigned char busnumber, unsigned char busaddress)
- array< int > [^] [GetPressure](#) (unsigned char busnumber, unsigned char busaddress, int n)
- int [GetPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRegulatorOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned char onoff)
- unsigned char [GetRegulatorOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)

- void [SetSollPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, int pressure)
- int [GetSollPressure](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRegulatorFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, int factor)
- int [GetRegulatorFactor](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- array< unsigned short > ^ [GetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress)
- int [GetPressureOffset](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- unsigned int [GetRegulatorStatus](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetRotatePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, short speed)
- short [GetRotatePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetMovePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned short speed, int position)
- void [GetMovePump](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% speed, [System::Runtime::InteropServices::Out]int% position)
- void [SetRegulationTimeouts](#) (unsigned char busnumber, unsigned char busaddress, unsigned short MaxSpeedWait, unsigned short MaxSignChange)
- void [GetRegulationTimeouts](#) (unsigned char busnumber, unsigned char busaddress, [System::Runtime::InteropServices::Out]unsigned short% MaxSpeedWait, [System::Runtime::InteropServices::Out]unsigned short% MaxSignChange)
- array< int > ^ [Get4ADC](#) (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ [Get4ADCAverage](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4DAC](#) (unsigned char busnumber, unsigned char busaddress, array< unsigned short > ^ dac)
- array< unsigned short > ^ [Get4DAC](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4ADCMode](#) (unsigned char busnumber, unsigned char busaddress, [PatchServAdcModeEnumNet](#) mode)
- [PatchServAdcModeEnumNet](#) [Get4ADCMode](#) (unsigned char busnumber, unsigned char busaddress)
- void [Set4ADCCatchampAverageShift](#) (unsigned char busnumber, unsigned char busaddress, unsigned int shift)
- unsigned int [Get4ADCCatchampAverageShift](#) (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ [Get2AnalogInput](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [Get2DigitalInput](#) (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ [GetADCs](#) (unsigned char busnumber, unsigned char busaddress, int n)
- array< unsigned short > ^ [GetADCsLoop](#) (unsigned char busnumber, unsigned char busaddress, int n)
- void [SetPiezoState](#) (unsigned char busnumber, unsigned char busaddress, int state)
- void [GetPiezoState](#) (unsigned char busnumber, unsigned char busaddress, [System::Runtime::InteropServices::Out]int% state, [System::Runtime::InteropServices::Out]int% reason)
- void [SetDACs](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index, array< unsigned short > ^ dac_times_voltages)
- array< unsigned short > ^ [GetDACs](#) (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void [SetSamplePeriode](#) (unsigned char busnumber, unsigned char busaddress, unsigned short periode)
- unsigned short [GetSamplePeriode](#) (unsigned char busnumber, unsigned char busaddress)
- void [StartSync](#) (unsigned char busnumber, unsigned char busaddress)
- unsigned short [GetSyncState](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacAmplitude](#) (unsigned char busnumber, unsigned char busaddress, unsigned short dacAmplitude)
- unsigned short [CatchAmpGetDacAmplitude](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacOffset](#) (unsigned char busnumber, unsigned char busaddress, short dacOffset)
- short [CatchAmpGetDacOffset](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcMean](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcValue](#) (unsigned char busnumber, unsigned char busaddress)

- int [CatchAmpGetAdcValueH](#) (unsigned char busnumber, unsigned char busaddress)
- int [CatchAmpGetAdcValueL](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetPwmEnable](#) (unsigned char busnumber, unsigned char busaddress, bool pwmEnable)
- bool [CatchAmpGetPwmEnable](#) (unsigned char busnumber, unsigned char busaddress)
- void [CatchAmpSetDacEnable](#) (unsigned char busnumber, unsigned char busaddress, bool dacEnable)
- bool [CatchAmpGetDacEnable](#) (unsigned char busnumber, unsigned char busaddress)
- int [TactSwitchGetState](#) (unsigned char busnumber, unsigned char busaddress)
- void [TactSwitchSetDisplay](#) (unsigned char busnumber, unsigned char busaddress, int Melody)

Additional Inherited Members

11.46.1 Constructor & Destructor Documentation

11.46.1.1 CMcsBus_SensorNet() [CMcsBus_SensorNet](#) (
[CMcsUsbNet](#)[^] *device*)

11.46.1.2 ~CMcsBus_SensorNet() [~CMcsBus_SensorNet](#) (
void)

11.46.2 Member Function Documentation

11.46.2.1 CatchAmpGetAdcMean() int [CatchAmpGetAdcMean](#) (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.46.2.2 CatchAmpGetAdcValue() int [CatchAmpGetAdcValue](#) (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.46.2.3 CatchAmpGetAdcValueH() int [CatchAmpGetAdcValueH](#) (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.46.2.4 CatchAmpGetAdcValueL() int CatchAmpGetAdcValueL (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.5 CatchAmpGetDacAmplitude() unsigned short CatchAmpGetDacAmplitude (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.6 CatchAmpGetDacEnable() bool CatchAmpGetDacEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.7 CatchAmpGetDacOffset() short CatchAmpGetDacOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.8 CatchAmpGetPwmEnable() bool CatchAmpGetPwmEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.9 CatchAmpSetDacAmplitude() void CatchAmpSetDacAmplitude (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *dacAmplitude*)

11.46.2.10 CatchAmpSetDacEnable() void CatchAmpSetDacEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 bool *dacEnable*)

11.46.2.11 CatchAmpSetDacOffset() void CatchAmpSetDacOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *dacOffset*)

11.46.2.12 CatchAmpSetPwmEnable() void CatchAmpSetPwmEnable (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 bool *pwmEnable*)

11.46.2.13 Get2AnalogInput() array<unsigned short> ^ Get2AnalogInput (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.14 Get2DigitalInput() unsigned short Get2DigitalInput (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.15 Get4ADC() array<int> ^ Get4ADC (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.16 Get4ADCAverage() array<int> ^ Get4ADCAverage (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.17 Get4ADCCatchampAverageShift() unsigned int Get4ADCCatchampAverageShift (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.18 Get4ADCMode() [PatchServAdcModeEnumNet](#) Get4ADCMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.19 Get4DAC() array<unsigned short> ^ Get4DAC (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.20 GetADCs() `array<unsigned short> ^ GetADCs (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `int n)`

11.46.2.21 GetADCsLoop() `array<unsigned short> ^ GetADCsLoop (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `int n)`

11.46.2.22 GetBubbleStatus() `unsigned short GetBubbleStatus (`
 `unsigned char busnumber,`
 `unsigned char busaddress)`

11.46.2.23 GetDACs() `array<unsigned short> ^ GetDACs (`
 `unsigned char busnumber,`
 `unsigned char busaddress,`
 `unsigned short index)`

11.46.2.24 GetDetectionThreshold() `unsigned short GetDetectionThreshold (`
 `unsigned char busnumber,`
 `unsigned char busaddress)`

11.46.2.25 GetDetectorValue() `unsigned short GetDetectorValue (`
 `unsigned char busnumber,`
 `unsigned char busaddress)`

11.46.2.26 GetLatency() `unsigned short GetLatency (`
 `unsigned char busnumber,`
 `unsigned char busaddress)`

11.46.2.27 GetLatencyCounter() `unsigned short GetLatencyCounter (`
 `unsigned char busnumber,`
 `unsigned char busaddress)`

11.46.2.28 GetMinimalThreshold() unsigned short GetMinimalThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.29 GetMovePump() void GetMovePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned short% *speed*,
 [System::Runtime::InteropServices::Out] int% *position*)

11.46.2.30 GetPiezoState() void GetPiezoState (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 [System::Runtime::InteropServices::Out] int% *state*,
 [System::Runtime::InteropServices::Out] int% *reason*)

11.46.2.31 GetPressure() [1/2] array<int> ^ GetPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *n*)

11.46.2.32 GetPressure() [2/2] int GetPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.33 GetPressureOffset() [1/2] array<unsigned short> ^ GetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.34 GetPressureOffset() [2/2] int GetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.35 GetRegulationTimeouts() void GetRegulationTimeouts (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 [System::Runtime::InteropServices::Out] unsigned short% *MaxSpeedWait*,
 [System::Runtime::InteropServices::Out] unsigned short% *MaxSignChange*)

11.46.2.36 GetRegulatorFactor() int GetRegulatorFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.37 GetRegulatorOnOff() unsigned char GetRegulatorOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.38 GetRegulatorStatus() unsigned int GetRegulatorStatus (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.39 GetRotatePump() short GetRotatePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.40 GetSamplePeriode() unsigned short GetSamplePeriode (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.41 GetSollPressure() int GetSollPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.42 GetSyncState() unsigned short GetSyncState (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.43 Set4ADCCatchampAverageShift() void Set4ADCCatchampAverageShift (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned int *shift*)

11.46.2.44 Set4ADCMode() void Set4ADCMode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 PatchServAdcModeEnumNet *mode*)

11.46.2.45 Set4DAC() void Set4DAC (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 array< unsigned short >^ *dac*)

11.46.2.46 SetDACs() void SetDACs (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 array< unsigned short >^ *dac_times_voltages*)

11.46.2.47 SetDetectionThreshold() void SetDetectionThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *threshold*)

11.46.2.48 SetLatency() void SetLatency (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *latency*)

11.46.2.49 SetMinimalThreshold() void SetMinimalThreshold (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *threshold*)

11.46.2.50 SetMovePump() void SetMovePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 unsigned short *speed*,
 int *position*)

11.46.2.51 SetPiezoState() void SetPiezoState (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *state*)

11.46.2.52 SetPressureOffset() void SetPressureOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*)

11.46.2.53 SetRegulationTimeouts() void SetRegulationTimeouts (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *MaxSpeedWait*,
 unsigned short *MaxSignChange*)

11.46.2.54 SetRegulatorFactor() void SetRegulatorFactor (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 int *factor*)

11.46.2.55 SetRegulatorOnOff() void SetRegulatorOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 unsigned char *onoff*)

11.46.2.56 SetRotatePump() void SetRotatePump (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 short *speed*)

11.46.2.57 SetSamplePeriode() void SetSamplePeriode (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *periode*)

11.46.2.58 SetSollPressure() void SetSollPressure (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *index*,
 int *pressure*)


11.46.2.59 StartSync() void StartSync (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.60 TactSwitchGetState() int TactSwitchGetState (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.46.2.61 TactSwitchSetDisplay() void TactSwitchSetDisplay (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *Melody*)

11.47 CMcsBus_TempSensorNet Class Reference

Inheritance diagram for CMcsBus_TempSensorNet:

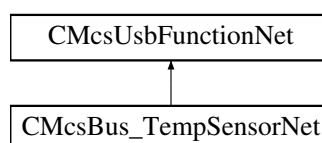

Public Member Functions

- [CMcsBus_TempSensorNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_TempSensorNet](#) (void)
- short [GetTemperatur](#) (unsigned char busnumber, unsigned char busaddress)
- short [GetTemperatur](#) (unsigned char busnumber, unsigned char busaddress, short index)
- void [SetNanoVoltsPerKelvin](#) (unsigned char busnumber, unsigned char busaddress, int nanovoltsperkelvin)
- int [GetNanoVoltsPerKelvin](#) (unsigned char busnumber, unsigned char busaddress)
- short [GetThermoVoltage](#) (unsigned char busnumber, unsigned char busaddress, short index)
- short [GetThermoTemp](#) (unsigned char busnumber, unsigned char busaddress, short index)
- void [SetThermoOffset](#) (unsigned char busnumber, unsigned char busaddress, short index, short offset)
- short [GetThermoOffset](#) (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members

11.47.1 Constructor & Destructor Documentation

11.47.1.1 [CMcsBus_TempSensorNet\(\)](#) [CMcsBus_TempSensorNet](#) (
 [CMcsUsbNet](#)[^] device)

11.47.1.2 [~CMcsBus_TempSensorNet\(\)](#) [~CMcsBus_TempSensorNet](#) (
 void)

11.47.2 Member Function Documentation

11.47.2.1 [GetNanoVoltsPerKelvin\(\)](#) int [GetNanoVoltsPerKelvin](#) (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.47.2.2 [GetTemperatur\(\)](#) [1/2] short [GetTemperatur](#) (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.47.2.3 [GetTemperatur\(\)](#) [2/2] short [GetTemperatur](#) (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *index*)

11.47.2.4 GetThermoOffset() short GetThermoOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *index*)

11.47.2.5 GetThermoTemp() short GetThermoTemp (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *index*)


11.47.2.6 GetThermoVoltage() short GetThermoVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *index*)

11.47.2.7 SetNanoVoltsPerKelvin() void SetNanoVoltsPerKelvin (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 int *nanovoltsperkelvin*)

11.47.2.8 SetThermoOffset() void SetThermoOffset (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 short *index*,
 short *offset*)

11.48 CMcsBus_VoltageModeNet Class Reference

Inheritance diagram for CMcsBus_VoltageModeNet:

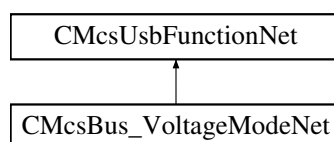

Public Member Functions

- [CMcsBus_VoltageModeNet](#) ([CMcsUsbNet](#)[^] device)
- [~CMcsBus_VoltageModeNet](#) (void)
- void [SetVMMMaxPositiveCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxPositiveCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxPositiveCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxNegativeCurrentEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short [GetVMMMaxNegativeCurrent](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxPositiveVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxPositiveVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxPositiveVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxNegativeVoltageEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMMMaxNegativeVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMMMaxNegativeVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMOOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, unsigned short status)
- unsigned short [GetVMOOutputOnOff](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void [SetVMVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short [GetVMVoltage](#) (unsigned char busnumber, unsigned char busaddress, unsigned char channel)

Additional Inherited Members

11.48.1 Constructor & Destructor Documentation

11.48.1.1 [CMcsBus_VoltageModeNet\(\)](#) [CMcsBus_VoltageModeNet](#) ([CMcsUsbNet](#)[^] device)

11.48.1.2 `~CMcsBus_VoltageModeNet()` `~CMcsBus_VoltageModeNet` (
 void)

11.48.2 Member Function Documentation

11.48.2.1 `GetVMMaxNegativeCurrent()` short `GetVMMaxNegativeCurrent` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.2 `GetVMMaxNegativeCurrentEeprom()` short `GetVMMaxNegativeCurrentEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.3 `GetVMMaxNegativeVoltage()` short `GetVMMaxNegativeVoltage` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.4 `GetVMMaxNegativeVoltageEeprom()` short `GetVMMaxNegativeVoltageEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.5 `GetVMMaxPositiveCurrent()` short `GetVMMaxPositiveCurrent` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.6 `GetVMMaxPositiveCurrentEeprom()` short `GetVMMaxPositiveCurrentEeprom` (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.7 GetVMMaxPositiveVoltage() short GetVMMaxPositiveVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.8 GetVMMaxPositiveVoltageEeprom() short GetVMMaxPositiveVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.9 GetVMOutputOnOff() unsigned short GetVMOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.10 GetVMVoltage() short GetVMVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*)

11.48.2.11 SetVMMaxNegativeCurrent() void SetVMMaxNegativeCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

11.48.2.12 SetVMMaxNegativeCurrentEeprom() void SetVMMaxNegativeCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

11.48.2.13 SetVMMaxNegativeVoltage() void SetVMMaxNegativeVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

11.48.2.14 SetVMMaxNegativeVoltageEeprom() void SetVMMaxNegativeVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

11.48.2.15 SetVMMaxPositiveCurrent() void SetVMMaxPositiveCurrent (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

11.48.2.16 SetVMMaxPositiveCurrentEeprom() void SetVMMaxPositiveCurrentEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *current*)

11.48.2.17 SetVMMaxPositiveVoltage() void SetVMMaxPositiveVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)


11.48.2.18 SetVMMaxPositiveVoltageEeprom() void SetVMMaxPositiveVoltageEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

11.48.2.19 SetVMOutputOnOff() void SetVMOutputOnOff (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 unsigned short *status*)

11.48.2.20 SetVMVoltage() void SetVMVoltage (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *channel*,
 short *voltage*)

11.49 CMcsBusNet Class Reference

Inheritance diagram for CMcsBusNet:

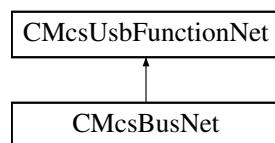

Public Member Functions

- [CMcsBusNet](#) ([CMcsUsbNet](#)[^] device)
- virtual [~CMcsBusNet](#) (void)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, short value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned int value)
- void [SetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, int value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned short% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]short% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned int% value)
- void [GetCommand](#) (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]int% value)
- void [SetBusAddressEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short [GetBusAddressEeprom](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetBusAddress](#) (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short [GetBusAddress](#) (unsigned char busnumber, unsigned char busaddress)
- void [CMcsBusNet::SetModeEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short [CMcsBusNet::GetModeEeprom](#) (unsigned char busnumber, unsigned char busaddress)
- void [CMcsBusNet::SetMode](#) (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short [CMcsBusNet::GetMode](#) (unsigned char busnumber, unsigned char busaddress)
- void [SetHWRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress, unsigned short revision)
- unsigned short [GetHWRevisionEeprom](#) (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.49.1 Constructor & Destructor Documentation

11.49.1.1 CMcsBusNet() [CMcsBusNet](#) (
[CMcsUsbNet](#)[^] device)

11.49.1.2 `~CMcsBusNet()` virtual `~CMcsBusNet` (
void) [virtual]

11.49.2 Member Function Documentation

11.49.2.1 `CMcsBusNet::GetMode()` unsigned short `CMcsBusNet::GetMode` (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.49.2.2 `CMcsBusNet::GetModeEeprom()` unsigned short `CMcsBusNet::GetModeEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.49.2.3 `CMcsBusNet::SetMode()` void `CMcsBusNet::SetMode` (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned short *mode*)

11.49.2.4 `CMcsBusNet::SetModeEeprom()` void `CMcsBusNet::SetModeEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*,
unsigned short *mode*)

11.49.2.5 `GetBusAddress()` unsigned short `GetBusAddress` (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.49.2.6 `GetBusAddressEeprom()` unsigned short `GetBusAddressEeprom` (
unsigned char *busnumber*,
unsigned char *busaddress*)

11.49.2.7 GetCommand() [1/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] int% *value*)

11.49.2.8 GetCommand() [2/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] short% *value*)

11.49.2.9 GetCommand() [3/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned int% *value*)

11.49.2.10 GetCommand() [4/4] void GetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 [System::Runtime::InteropServices::Out] unsigned short% *value*)

11.49.2.11 GetHWRevisionEeprom() unsigned short GetHWRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*)

11.49.2.12 SetBusAddress() void SetBusAddress (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *newaddress*)

11.49.2.13 SetBusAddressEeprom() void SetBusAddressEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *newaddress*)

11.49.2.14 SetCommand() [1/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 int *value*)

11.49.2.15 SetCommand() [2/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 short *value*)

11.49.2.16 SetCommand() [3/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned int *value*)


11.49.2.17 SetCommand() [4/4] void SetCommand (
 unsigned char *command*,
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned char *axis*,
 unsigned short *value*)

11.49.2.18 SetHWRevisionEeprom() void SetHWRevisionEeprom (
 unsigned char *busnumber*,
 unsigned char *busaddress*,
 unsigned short *revision*)

11.50 CMcsUsbDacqNet Class Reference

Base class for data acquisition devices.

Inheritance diagram for CMcsUsbDacqNet:

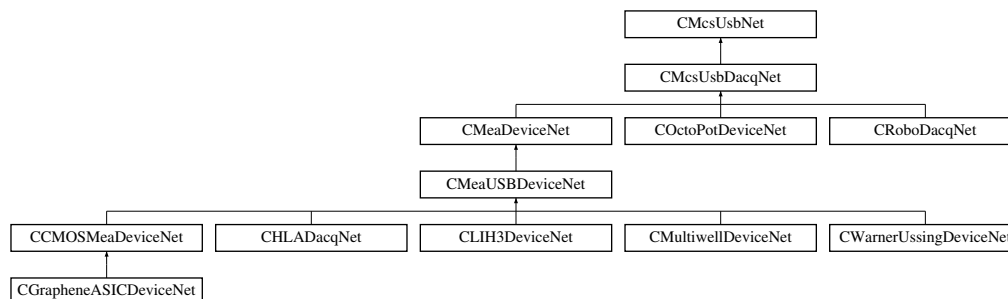

Classes

- class [CHWInfo](#)
Class to provide hardware information about the device.

Public Member Functions

- [CMcsUsbDacqNet](#) ()
- [~CMcsUsbDacqNet](#) ()
- [uint32_t GetErrorMessage](#) ([System::Runtime::InteropServices::Out]String^% errorString, [System::Runtime::InteropServices::Out]int% info)
- [virtual uint32_t GetVoltageRangeIndex](#) (unsigned int virtualDevice)
- [virtual void SetVoltageRangeByIndex](#) (int32_t voltageRangeIndex, unsigned int virtualDevice)
Sets the voltage range on devices which support multiple voltage ranges.
- [virtual void SetVoltageRangeInMicroVolt](#) (int32_t voltageRange, unsigned int virtualDevice)
Sets the voltage range on devices which support multiple voltage ranges.
- [virtual int32_t GetVoltageRangeInMicroVolt](#) (unsigned int virtualDevice)
Gets the currently selected voltage range on devices which support multiple voltage ranges.
- [virtual int32_t GetVoltageRangeInMilliVolt](#) ()
Gets the currently selected voltage range on devices which support multiple voltage ranges.
- [virtual void SetDataMode](#) (DataModeEnumNet dataMode, unsigned int virtualDevice)
Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.
- [virtual DataModeEnumNet GetDataMode](#) (unsigned int virtualDevice)
Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.
- [void SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- [void SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, W2100DigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- [void SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, SCUDigitalSourceEnumNet source, int bitnumber_offset)
Sets the function/source of an digital output bit.
- [void SetDigitalSource](#) (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MEA2100_256DigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

- template<typename digitalsourceenum >
void **SetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSource< digitalsourceenum >^ source, int bitnumber_offset)

Sets the function/source of an digital output bit.

- void **GetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

- void **GetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

- void **GetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]SCUDigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

- void **GetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]MEA2100_256DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

- template<typename digitalsourceenum >
void **GetDigitalSource** (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::InteropServices::Out]DigitalSource< digitalsourceenum >^% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

- virtual AdapterTypeEnumNet **GetAdapterType** ()

Gets the adapter which is connected to the MEA2100 device.

- virtual MeaLayoutEnumNet **GetMeaLayout** ()

Gets the MEA layout which is connected to the MEA2100 device.

- virtual uint32_t **GetAdcDataFormat** (uint32_t virtualDevice)

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

- virtual uint32_t **GetAnalogValueUnit** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] AnalogUnitEnumNet% unit)
- virtual uint32_t **GetResolutionPerDigit** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% res, [System::Runtime::InteropServices::Out] int% resUnit)
- virtual uint32_t **GetAdcZero** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% adcz)
- virtual uint32_t **GetHardwareMinRange** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t **GetHardwareMaxRange** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t **GetDataFormat** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)
- virtual uint32_t **GetNumberOfDataBits** (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)

Get the real number of data bits.

- virtual void **SetSamplerate** (int32_t rate, unsigned int oversample, unsigned int virtualDevice)

Sets the sampling frequency of the device.

- virtual int32_t **GetSamplerate** (unsigned int virtualDevice)

Gets the sampling frequency of the device.

- virtual uint32_t **GetMaxSamplingFrequency** (int virtualDevice)

Gets the maximal sampling frequency of the device.

- virtual uint32_t **GetMinSamplingFrequencyStepsize** ()

Gets the minimal sampling frequency step size increment value of the device.

- virtual int32_t [GetChannelsInBlock](#) (unsigned int virtualDevice)

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

- virtual void [GetChannelLayout](#) ([System::Runtime::InteropServices::Out]int% AnalogChannels, [System::Runtime::InteropServices::Out]int% DigitalChannels, [System::Runtime::InteropServices::Out]int% ChecksumChannels, [System::Runtime::InteropServices::Out]int% TimestampChannels, [System::Runtime::InteropServices::Out]int% ChannelsInBlock, unsigned int virtualDevice)

- virtual void [SendStartDacq](#) ()

Start sampling.

- virtual void [SendStartDacq](#) (int VirtualDacqMap)

Start sampling.

- virtual void [SendStartStgAndDacq](#) (uint32_t trigger_map, int VirtualDacqMap)

Start sampling together with the STG.

- virtual void [SendStopDacq](#) ()

Stop sampling.

- virtual void [SendStopDacq](#) (int VirtualDacqMap)

Stop sampling.

Parameters

VirtualDacqMap	
----------------	--

- virtual void [SendStopStgAndDacq](#) (uint32_t trigger_map, int VirtualDacqMap)

Stop sampling together with the STG.

- virtual void [SendStopStgAndDacqWithOptions](#) (uint32_t trigger_map, int VirtualDacqMap, int options)

Stop sampling together with the STG and options.

- virtual void [StartLoop](#) ()

Start the data acquisition thread.

- virtual void [StartLoop](#) (int32_t timeout)

Start the data acquisition thread.

- virtual void [StartLoop](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb)

Start the data acquisition thread.

- virtual void [StartLoop](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread.

- virtual void [StopLoop](#) ()

- virtual void [ClearBuffers](#) ()

- virtual void [StartDacq](#) ()

Start the data acquisition thread and sampling.

- virtual void [StartDacq](#) (int32_t timeout)

Start the data acquisition thread and sampling.

- virtual void [StartDacq](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb)

Start the data acquisition thread and sampling.

- virtual void [StartDacq](#) (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread and sampling.

- virtual void [StopDacq](#) ()

Stop the data acquisition thread and sampling.

- virtual void [StopDacq](#) (uint32_t virtualDevice)

Stop the data acquisition thread and sampling.

- virtual uint32_t [SetPoti](#) (uint32_t channel, uint32_t value, bool write_nvram)
- virtual uint32_t [GetPoti](#) (uint32_t channel, [System::Runtime::InteropServices::Out]uint32_t% value)
- virtual [CFilterPropertyNet](#) ^ [GetFilterProperty](#) ([DacqGroupChannelEnumNet](#) GroupID, unsigned int index)
- virtual array< [CFilterPropertyNet](#) ^ > ^ [CMcsUsbDacqNet::GetFilterProperties](#) ([DacqGroupChannelEnumNet](#) GroupID)
- int [GetChannelDataFillSize](#) ()
- virtual void [SetSelectedChannels](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void [SetSelectedChannels](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedChannels](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void [SetSelectedChannels](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedData](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void [SetSelectedData](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedData](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void [SetSelectedData](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize)
- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual int [AddSelectedChannelsQueue](#) (int nByteOffset, int nChannelOffset, array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize)
- virtual void [SetSelectedChannelsQueue](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual void [SetSelectedChannelsQueue](#) (int nChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual void [SetSelectedChannelsQueue](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_ReadFramesDict... with handle = 0 to read the data.

- virtual void [SetSelectedChannelsQueue](#) (array< bool > ^ selectedChannels, int queuesize, int threshold, [SampleSizeNet](#) samplesize, [SampleDstSizeNet](#) sampleDstSize, int ChannelsInBlock)
- virtual uint32_t [ChannelBlock_AvailFrames](#) (int handle)

Get the number of sample frames already available in the FIFO.

- virtual uint32_t [ChannelBlock_AvailFrames](#) (int handle, int queue)

- virtual array< uint16_t > ^ [ChannelBlock_ReadFramesUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint16_t data format
- virtual void [ChannelBlock_ReadFramesUI16](#) (int handle, array< uint16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint16_t data format
- virtual array< int16_t > ^ [ChannelBlock_ReadFramesI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int16_t data format
- virtual void [ChannelBlock_ReadFramesI16](#) (int handle, array< int16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int16_t data format
- virtual array< uint32_t > ^ [ChannelBlock_ReadFramesUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format
- virtual void [ChannelBlock_ReadFramesUI32](#) (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format
- virtual array< int32_t > ^ [ChannelBlock_ReadFramesI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format
- virtual void [ChannelBlock_ReadFramesI32](#) (int handle, array< int32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format
- virtual array< array< uint16_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< uint16_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI16](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< int16_t >^> ^ [ChannelBlock_ReadAsFrameArrayI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< int16_t >^> ^ [ChannelBlock_ReadAsFrameArrayI16](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< uint32_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< uint32_t >^> ^ [ChannelBlock_ReadAsFrameArrayUI32](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< int32_t >^> ^ [ChannelBlock_ReadAsFrameArrayI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual array< array< int32_t >^> ^ [ChannelBlock_ReadAsFrameArrayI32](#) (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue as array of uint16_t data frame arrays
- virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ [ChannelBlock_ReadFramesDictUI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

- virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ [ChannelBlock_ReadFramesDictI16](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ [ChannelBlock_ReadFramesDictUI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ [ChannelBlock_ReadFramesDictI32](#) (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ [GetGroupChannelDataUI16](#) ([DacqGroupChannelEnumNet](#) group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ [GetGroupChannelDataI16](#) ([DacqGroupChannelEnumNet](#) group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ [GetGroupChannelDataUI32](#) ([DacqGroupChannelEnumNet](#) group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ [GetGroupChannelDataI32](#) ([DacqGroupChannelEnumNet](#) group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number
- void [SetupGroupDacqQueue](#) (int queuesize, int threshold)
- void [SetupGroupDacqQueue](#) (int queuesize, int threshold, unsigned int virtualDevice)
- [CHWInfo](#) ^ [HWInfo](#) ()

Static Public Attributes

- static const int [Error_Callback_Queue_Full](#) = 0x100
- static const int [Error_Callback_Aquisition_Stopped](#) = 0x200
- static const int [Error_Callback_Packet_Error](#) = 1
- static const int [Error_Callback_RingQueue_Full](#) = 3
- static const int [Error_Callback_Frames_Lost](#) = 4
- static const int [Error_Callback_Data_lost](#) = 5

Properties

- virtual int [Samplerate](#) [get, set]
The sampling frequency of the device in Hz.

Events

- [OnChannelData](#) ^ [ChannelDataEvent](#) [add, remove, raise]
- [OnError](#) ^ [ErrorEvent](#) [add, remove, raise]

Additional Inherited Members

11.50.1 Detailed Description

Base class for data acquisition devices.

11.50.2 Constructor & Destructor Documentation

11.50.2.1 CMcsUsbDacqNet() `CMcsUsbDacqNet ()`

11.50.2.2 ~CMcsUsbDacqNet() `~CMcsUsbDacqNet ()`

11.50.3 Member Function Documentation

11.50.3.1 AddSelectedChannelsQueue() [1/4] `virtual int AddSelectedChannelsQueue (`
`int nByteOffset,`
`int nChannelOffset,`
`array< bool >^ selectedChannels,`
`int queuesize,`
`int threshold,`
`SampleSizeNet samplesize) [virtual]`

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read↔FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>nByteOffset</i>	Number of bytes to start with.
--------------------	--------------------------------

Parameters

<i>nChannelOffset</i>	Number of channel to start with (counted in samplesize bytes).
-----------------------	--

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Returns

The handle to the Queue.

11.50.3.2 AddSelectedChannelsQueue() [2/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*) [virtual]

11.50.3.3 AddSelectedChannelsQueue() [3/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*) [virtual]

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use `ChannelBlock_Read↔FramesDict...` with `handle = 0` to read the data.

When using 32 bit data format, `ChannelsInBlock` is still the number of 16 bit channels per frame, as obtained from `GetChannelsInBlock`, while `nChannels` is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format $nChannels = ChannelsInBlock/2$

Parameters

<i>nByteOffset</i>	Number of bytes to start with.
--------------------	--------------------------------

Parameters

<i>nChannelOffset</i>	Number of channel to start with (counted in samplesize bytes).
-----------------------	--

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Returns

The handle to the Queue.

11.50.3.4 AddSelectedChannelsQueue() [4/4] virtual int AddSelectedChannelsQueue (
 int *nByteOffset*,
 int *nChannelOffset*,
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*) [virtual]

11.50.3.5 ChannelBlock_AvailFrames() [1/2] virtual uint32_t ChannelBlock_AvailFrames (
 int *handle*) [virtual]

Get the number of sample frames already available in the FIFO.

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Returns

Number of sample frames available in the FIFO.

11.50.3.6 ChannelBlock_AvailFrames() [2/2] virtual uint32_t ChannelBlock_AvailFrames (
 int *handle*,
 int *queue*) [virtual]

11.50.3.7 ChannelBlock_ReadAsFrameArrayI16() [1/2] virtual array<array<int16_t>> ^ Channel↔
 Block_ReadAsFrameArrayI16 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int16_t frame arrays.

```
11.50.3.8 ChannelBlock_ReadAsFrameArrayI16() [2/2]  virtual array<array<int16_t>^> ^ Channel↔
Block_ReadAsFrameArrayI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int16_t frame arrays.

11.50.3.9 ChannelBlock_ReadAsFrameArrayI32() [1/2] `virtual array<array<int32_t>^> ^ Channel↔
Block_ReadAsFrameArrayI32 (`
`int handle,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int32_t frame arrays.

11.50.3.10 ChannelBlock_ReadAsFrameArrayI32() [2/2] `virtual array<array<int32_t>^> ^ ChannelBlock_ReadAsFrameArrayI32 (`
`int handle,`
`int queue,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of int32_t frame arrays.

11.50.3.11 ChannelBlock_ReadAsFrameArrayUI16() [1/2] `virtual array<array<uint16_t>^> ^ ChannelBlock_ReadAsFrameArrayUI16 (`
`int handle,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint16_t frame arrays.

11.50.3.12 ChannelBlock_ReadAsFrameArrayUI16() [2/2] virtual array<array<uint16_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI16 (  
    int handle,  
    int queue,  
    int frames,  
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint16_t frame arrays.

11.50.3.13 ChannelBlock_ReadAsFrameArrayUI32() [1/2] `virtual array<array<uint32_t>^> ^`
`ChannelBlock_ReadAsFrameArrayUI32 (`
 `int handle,`
 `int frames,`
 `[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint32_t frame arrays.

11.50.3.14 ChannelBlock_ReadAsFrameArrayUI32() [2/2] `virtual array<array<uint32_t>^> ^`
`ChannelBlock_ReadAsFrameArrayUI32 (`
 `int handle,`
 `int queue,`
 `int frames,`
 `[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
---------------	--

Parameters

<i>queue</i>	Number of the sub queue.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of uint32_t frame arrays.

11.50.3.15 ChannelBlock_ReadFramesDictI16() virtual System::Collections::Generic::Dictionary<int, array<int16_t>^> ^ ChannelBlock_ReadFramesDictI16 (int *handle*, int *frames*, [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int16_t arrays and hardware channel as key.

```
11.50.3.16 ChannelBlock_ReadFramesDictI32() virtual System::Collections::Generic::Dictionary<int,
array<int32_t>^> ^ ChannelBlock_ReadFramesDictI32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int32_t arrays and hardware channel as key.

```
11.50.3.17 ChannelBlock_ReadFramesDictUI16() virtual System::Collections::Generic::Dictionary<int,
array<uint16_t>^> ^ ChannelBlock_ReadFramesDictUI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint16_t arrays and hardware channel as key.

11.50.3.18 ChannelBlock_ReadFramesDictUI32() `virtual System::Collections::Generic::Dictionary<int, array<uint32_t>^> ^ ChannelBlock_ReadFramesDictUI32 (int handle, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>handle</i>	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint32_t arrays and hardware channel as key.

11.50.3.19 ChannelBlock_ReadFramesI16() `[1/2] virtual void ChannelBlock_ReadFramesI16 (int handle, array< int16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in int16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

```
11.50.3.20 ChannelBlock_ReadFramesI16() [2/2] virtual array<int16_t> ^ ChannelBlock_Read↵
FramesI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

```
11.50.3.21 ChannelBlock_ReadFramesI32() [1/2] virtual void ChannelBlock_ReadFramesI32 (
    int handle,
    array< int32_t >^ buffer,
    int frames_pos,
```



```
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

11.50.3.22 ChannelBlock_ReadFramesI32() [2/2] virtual array<int32_t> ^ ChannelBlock_ReadFramesI32 (

```
int handle,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

11.50.3.23 ChannelBlock_ReadFramesUI16() [1/2] `virtual void ChannelBlock_ReadFramesUI16 (`
`int handle,`
`array< uint16_t >^ buffer,`
`int frames_pos,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

11.50.3.24 ChannelBlock_ReadFramesUI16() [2/2] `virtual array<uint16_t> ^ ChannelBlock_ReadFramesUI16 (`
`int handle,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

Returns

Array of data from the device.

11.50.3.25 ChannelBlock_ReadFramesUI32() [1/2] `virtual void ChannelBlock_ReadFramesUI32 (`
`int handle,`
`array< uint32_t >^ buffer,`
`int frames_pos,`
`int frames,`
`[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>buffer</i>	Buffer to put the data from the device in.
<i>frames_pos</i>	Position in buffer where to put the data.
<i>frames</i>	Number of sample frames to read.

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

11.50.3.26 ChannelBlock_ReadFramesUI32() [2/2] virtual array<uint32_t> ^ ChannelBlock_Read↔
FramesUI32 (
 int *handle*,
 int *frames*,
 [System::Runtime::InteropServices::Out] int % *frames_ret*) [virtual]

Read data from a FIFO queue in uint32_t data format

Parameters

<i>handle</i>	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
---------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
---------------	----------------------------------

Parameters

<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.
-------------------	--

11.50.3.27 ClearBuffers() virtual void ClearBuffers () [virtual]

11.50.3.28 CMcsUsbDacqNet::GetFilterProperties() virtual array<CFilterPropertyNet^> ^ CMcs↔
UsbDacqNet::GetFilterProperties (
 DacqGroupChannelEnumNet *GroupID*) [virtual]

11.50.3.29 GetAdapterType() virtual AdapterTypeEnumNet GetAdapterType () [virtual]

Gets the adapter which is connected to the MEA2100 device.

Returns

AdapterTypeEnumNet which enumerates the possible adapters.

11.50.3.30 GetAdcDataFormat() virtual uint32_t GetAdcDataFormat (
 uint32_t virtualDevice) [virtual]

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

Returns

The data format in bits.

11.50.3.31 GetAdcZero() virtual uint32_t GetAdcZero (
 uint32_t virtualDevice,
 DacqGroupChannelEnumNet group,
 [System::Runtime::InteropServices::Out] int% adcz) [virtual]

11.50.3.32 GetAnalogValueUnit() virtual uint32_t GetAnalogValueUnit (
 uint32_t virtualDevice,
 DacqGroupChannelEnumNet group,
 [System::Runtime::InteropServices::Out] AnalogUnitEnumNet% unit) [virtual]

11.50.3.33 GetChannelDataFillSize() int GetChannelDataFillSize ()

11.50.3.34 GetChannelLayout() virtual void GetChannelLayout (
 [System::Runtime::InteropServices::Out] int% AnalogChannels,
 [System::Runtime::InteropServices::Out] int% DigitalChannels,
 [System::Runtime::InteropServices::Out] int% ChecksumChannels,
 [System::Runtime::InteropServices::Out] int% TimestampChannels,
 [System::Runtime::InteropServices::Out] int% ChannelsInBlock,
 unsigned int virtualDevice) [virtual]

11.50.3.35 GetChannelsInBlock() virtual int32_t GetChannelsInBlock (
 unsigned int virtualDevice) [virtual]

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

Returns

Number of 16 bit datawords per sample frame.

11.50.3.36 GetDataFormat() virtual uint32_t GetDataFormat (
 uint32_t virtualDevice,
 DacqGroupChannelEnumNet group,
 [System::Runtime::InteropServices::Out] int% numberOfBits) [virtual]

11.50.3.37 GetDataMode() virtual DataModeEnumNet GetDataMode (
 unsigned int virtualDevice) [virtual]

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------

Returns

DataModeEnumNet which enumerates the possible data modes.

11.50.3.38 GetDigitalSource() [1/5] void GetDigitalSource (
 DigitalTargetEnumNet digitaltarget,
 int32_t NrChannel,
 [System::Runtime::InteropServices::Out] DigitalSource< digitalsourceenum >^%
 source,
 [System::Runtime::InteropServices::Out] int% bitnumber_offset)

Gets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.39 GetDigitalSource() [2/5] void GetDigitalSource (
 DigitalTargetEnumNet digitaltarget,
 int32_t NrChannel,
 [System::Runtime::InteropServices::Out] DigitalSourceEnumNet% source,
 [System::Runtime::InteropServices::Out] int% bitnumber_offset)

Gets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.40 GetDigitalSource() [3/5] `void GetDigitalSource (
 DigitalTargetEnumNet digitaltarget,
 int32_t NrChannel,
 [System::Runtime::InteropServices::Out] MEA2100_256DigitalSourceEnumNet% source,
 [System::Runtime::InteropServices::Out] int% bitnumber_offset)`

Gets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.41 GetDigitalSource() [4/5] `void GetDigitalSource (
 DigitalTargetEnumNet digitaltarget,
 int32_t NrChannel,
 [System::Runtime::InteropServices::Out] SCUDigitalSourceEnumNet% source,
 [System::Runtime::InteropServices::Out] int% bitnumber_offset)`

Gets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.42 GetDigitalSource() [5/5] `void GetDigitalSource (
 DigitalTargetEnumNet digitaltarget,
 int32_t NrChannel,
 [System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% source,
 [System::Runtime::InteropServices::Out] int% bitnumber_offset)`

Gets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

<i>digitaltarget</i>	The digital target to query.
<i>NrChannel</i>	The channel/bit of target to query.
<i>source</i>	The source/function assignd to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.43 GetErrorMessage() uint32_t GetErrorMessage (

```
[System::Runtime::InteropServices::Out] String^% errorString,
[System::Runtime::InteropServices::Out] int% info )
```

11.50.3.44 GetFilterProperty() virtual CFilterPropertyNet ^ GetFilterProperty (

```
DacqGroupChannelEnumNet GroupID,
unsigned int index ) [virtual]
```

11.50.3.45 GetGroupChannelDataI16() virtual System::Collections::Generic::Dictionary<int, array<int16_t>^> ^ GetGroupChannelDataI16 (

```
DacqGroupChannelEnumNet group,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int16_t arrays and hardware channel as key.

11.50.3.46 GetGroupChannelDataI32() virtual System::Collections::Generic::Dictionary<int, array<int32_t>^> ^ GetGroupChannelDataI32 (

```
DacqGroupChannelEnumNet group,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of int32_t arrays and hardware channel as key.

11.50.3.47 GetGroupChannelDataUI16() `virtual System::Collections::Generic::Dictionary<int, array<uint16_t>>^> ^ GetGroupChannelDataUI16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]`

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint16_t arrays and hardware channel as key.

11.50.3.48 GetGroupChannelDataUI32() virtual System::Collections::Generic::Dictionary<int, array<uint32_t>^> ^ GetGroupChannelDataUI32 (
 DacqGroupChannelEnumNet group,
 int frames,
 [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

<i>group</i>	Group selector supported by the device.
--------------	---

Parameters

<i>frames</i>	Number of sample frames to read.
<i>frames_ret</i>	Number of sample frames which were read, might be smaller than frames.

Returns

Dictionary of uint32_t arrays and hardware channel as key.

11.50.3.49 GetHardwareMaxRange() virtual uint32_t GetHardwareMaxRange (
 uint32_t virtualDevice,
 DacqGroupChannelEnumNet group,
 [System::Runtime::InteropServices::Out] int% r,
 [System::Runtime::InteropServices::Out] int% rUnit) [virtual]

11.50.3.50 GetHardwareMinRange() virtual uint32_t GetHardwareMinRange (
 uint32_t virtualDevice,
 DacqGroupChannelEnumNet group,
 [System::Runtime::InteropServices::Out] int% r,
 [System::Runtime::InteropServices::Out] int% rUnit) [virtual]

11.50.3.51 GetMaxSamplingFrequency() virtual uint32_t GetMaxSamplingFrequency (
 int virtualDevice) [virtual]

Gets the maximal sampling frequency of the device.

Returns

Sampling frequency in Hz.

11.50.3.52 GetMeaLayout() virtual [MeaLayoutEnumNet](#) GetMeaLayout () [virtual]

Gets the MEA layout which is connected to the MEA2100 device.

Returns

MeaLayoutEnumNet which enumerates the MEA types.

11.50.3.53 GetMinSamplingFrequencyStepsize() virtual uint32_t GetMinSamplingFrequencyStepsize () [virtual]

Gets the minimal sampling frequency step size increment value of the device.

Returns

Sampling frequency step size in Hz.

11.50.3.54 GetNumberOfDataBits() virtual uint32_t GetNumberOfDataBits (uint32_t *virtualDevice*, [DacqGroupChannelEnumNet](#) *group*, [System::Runtime::InteropServices::Out] int% *numberOfBits*) [virtual]

Get the real number of data bits.

This value may be different from the value returned by GetDataFormat, e.g. in MC_Card the data are shifted 2 bits so the real number is 14 while the data format is 16 bits

11.50.3.55 GetPoti() virtual uint32_t GetPoti (uint32_t *channel*, [System::Runtime::InteropServices::Out] uint32_t% *value*) [virtual]

11.50.3.56 GetResolutionPerDigit() virtual uint32_t GetResolutionPerDigit (uint32_t *virtualDevice*, [DacqGroupChannelEnumNet](#) *group*, [System::Runtime::InteropServices::Out] int% *res*, [System::Runtime::InteropServices::Out] int% *resUnit*) [virtual]

11.50.3.57 GetSamplerate() virtual int32_t GetSamplerate (unsigned int *virtualDevice*) [virtual]

Gets the sampling frequency of the device.

Returns

Sampling frequency in Hz.

11.50.3.58 GetVoltageRangeIndex() `virtual uint32_t GetVoltageRangeIndex (unsigned int virtualDevice) [virtual]`

11.50.3.59 GetVoltageRangeInMicroVolt() `virtual int32_t GetVoltageRangeInMicroVolt (unsigned int virtualDevice) [virtual]`

Gets the currently selected voltage range on devices which support multiple voltage ranges.

Returns

The Voltage Range in uV.

11.50.3.60 GetVoltageRangeInMilliVolt() `virtual int32_t GetVoltageRangeInMilliVolt () [virtual]`

Gets the currently selected voltage range on devices which support multiple voltage ranges.

Returns

The rounded Voltage Range in mV.

11.50.3.61 HWInfo() `CHWInfo ^ HWInfo ()`

11.50.3.62 SendStartDacq() [1/2] `virtual void SendStartDacq () [virtual]`

Start sampling.

11.50.3.63 SendStartDacq() [2/2] `virtual void SendStartDacq (int VirtualDacqMap) [virtual]`

Start sampling.

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

11.50.3.64 SendStartStgAndDacq() virtual void SendStartStgAndDacq (
 uint32_t *trigger_map*,
 int *VirtualDacqMap*) [virtual]

Start sampling together with the STG.

Parameters

<i>trigger_map</i>	
<i>VirtualDacqMap</i>	

11.50.3.65 SendStopDacq() [1/2] virtual void SendStopDacq () [virtual]

Stop sampling.

11.50.3.66 SendStopDacq() [2/2] virtual void SendStopDacq (
 int *VirtualDacqMap*) [virtual]

Stop sampling.

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

11.50.3.67 SendStopStgAndDacq() virtual void SendStopStgAndDacq (
 uint32_t *trigger_map*,
 int *VirtualDacqMap*) [virtual]

Stop sampling together with the STG.

Parameters

<i>trigger_map</i>	
--------------------	--

11.50.3.68 SendStopStgAndDacqWithOptions() virtual void SendStopStgAndDacqWithOptions (
 uint32_t *trigger_map*,

```
int VirtualDacqMap,
int options ) [virtual]
```

Stop sampling together with the STG and options.

Parameters

<i>trigger_map</i>	
--------------------	--

Parameters

<i>options</i>	
----------------	--

Parameters

<i>VirtualDacqMap</i>	
-----------------------	--

11.50.3.69 SetDataMode() `virtual void SetDataMode (`
`DataModeEnumNet dataMode,`
`unsigned int virtualDevice) [virtual]`

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

<i>dataMode</i>	DataModeEnumNet enumerates the possible data modes.
<i>virtualDevice</i>	Virtual device to use.

11.50.3.70 SetDigitalSource() [1/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`DigitalSource< digitalsourceenum >^ source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.71 SetDigitalSource() [2/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`DigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.72 SetDigitalSource() [3/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`MEA2100_256DigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.73 SetDigitalSource() [4/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`

```

int32_t NrChannel,
SCUDigitalSourceEnumNet source,
int bitnumber_offset )

```

Sets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.74 SetDigitalSource() [5/5] `void SetDigitalSource (`
`DigitalTargetEnumNet digitaltarget,`
`int32_t NrChannel,`
`W2100DigitalSourceEnumNet source,`
`int bitnumber_offset)`

Sets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

<i>digitaltarget</i>	The digital target to change.
<i>NrChannel</i>	The channel/bit of target to change.
<i>source</i>	The source/function to assign to the digital target.
<i>bitnumber_offset</i>	An offset / bit number with the source/function.

11.50.3.75 SetPoti() `virtual uint32_t SetPoti (`
`uint32_t channel,`
`uint32_t value,`
`bool write_nvram) [virtual]`

11.50.3.76 SetSamplerate() `virtual void SetSamplerate (`
`int32_t rate,`
`unsigned int oversample,`
`unsigned int virtualDevice) [virtual]`

Sets the sampling frequency of the device.

Parameters

<i>rate</i>	Sampling frequency in Hz.
-------------	---------------------------

11.50.3.77 SetSelectedChannels() [1/4] virtual void SetSelectedChannels (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be divided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

11.50.3.78 SetSelectedChannels() [2/4] virtual void SetSelectedChannels (
array< bool >^ *selectedChannels*,
int *queuesize*,
int *threshold*,
SampleSizeNet *samplesize*,
SampleDstSizeNet *sampleDstSize*,
int *ChannelsInBlock*) [virtual]

11.50.3.79 SetSelectedChannels() [3/4] virtual void SetSelectedChannels (
int *nChannels*,
int *queuesize*,
int *threshold*,
SampleSizeNet *samplesize*,
int *ChannelsInBlock*) [virtual]

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be divided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of samples frames the FIFO must acquire before the callback function is called.
------------------	--

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

11.50.3.80 SetSelectedChannels() [4/4] virtual void SetSelectedChannels (

```

    int nChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    SampleDstSizeNet sampleDstSize,
    int ChannelsInBlock ) [virtual]

```

11.50.3.81 SetSelectedChannelsQueue() [1/4] virtual void SetSelectedChannelsQueue (

```

    array< bool >^ selectedChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    int ChannelsInBlock ) [virtual]

```

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_↔ ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

11.50.3.82 SetSelectedChannelsQueue() [2/4] virtual void SetSelectedChannelsQueue (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

11.50.3.83 SetSelectedChannelsQueue() [3/4] virtual void SetSelectedChannelsQueue (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_↔
ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
-------------------	--

Parameters

<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.
------------------------	---

11.50.3.84 SetSelectedChannelsQueue() [4/4] virtual void SetSelectedChannelsQueue (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

11.50.3.85 SetSelectedData() [1/4] virtual void SetSelectedData (
 array< bool >^ *selectedChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 int *ChannelsInBlock*) [virtual]

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>selectedChannels</i>	List of channels to be collected in the FIFO.
-------------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

11.50.3.86 SetSelectedData() [2/4] virtual void SetSelectedData (
 array< bool >^ selectedChannels,
 int queuesize,
 int threshold,
 SampleSizeNet samplesize,
 SampleDstSizeNet sampleDstSize,
 int ChannelsInBlock) [virtual]

11.50.3.87 SetSelectedData() [3/4] virtual void SetSelectedData (
 int nChannels,
 int queuesize,
 int threshold,
 SampleSizeNet samplesize,
 int ChannelsInBlock) [virtual]

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

<i>nChannels</i>	Number of channels to be collected in the FIFO.
------------------	---

Parameters

<i>queuesize</i>	Size of sample frames the FIFO can hold.
------------------	--

Parameters

<i>threshold</i>	Number of sample frames the FIFO must acquire before the callback function is called.
------------------	---

Parameters

<i>samplesize</i>	size of the datawords, either 16 or 32bit.
<i>ChannelsInBlock</i>	value obtained from GetChannelsInBlock.

11.50.3.88 SetSelectedData() [4/4] virtual void SetSelectedData (
 int *nChannels*,
 int *queuesize*,
 int *threshold*,
 SampleSizeNet *samplesize*,
 SampleDstSizeNet *sampleDstSize*,
 int *ChannelsInBlock*) [virtual]

11.50.3.89 SetupGroupDacqQueue() [1/2] void SetupGroupDacqQueue (
 int *queuesize*,
 int *threshold*)

11.50.3.90 SetupGroupDacqQueue() [2/2] void SetupGroupDacqQueue (
 int *queuesize*,
 int *threshold*,
 unsigned int *virtualDevice*)

11.50.3.91 SetVoltageRangeByIndex() virtual void SetVoltageRangeByIndex (
 int32_t *voltageRangeIndex*,
 unsigned int *virtualDevice*) [virtual]

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

<i>voltageRangeIndex</i>	Voltage Range to use as index, smaller values are larger voltage ranges.
--------------------------	--

11.50.3.92 SetVoltageRangeInMicroVolt() `virtual void SetVoltageRangeInMicroVolt (`
 `int32_t voltageRange,`
 `unsigned int virtualDevice) [virtual]`

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

<i>voltageRange</i>	Voltage Range to use in μ V.
---------------------	----------------------------------

This replaces SetVoltageRange, where the value of the range was in mV!

11.50.3.93 StartDacq() [1/4] `virtual void StartDacq () [virtual]`

Start the data acquisition thread and sampling.

11.50.3.94 StartDacq() [2/4] `virtual void StartDacq (`
 `int32_t timeout) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

11.50.3.95 StartDacq() [3/4] `virtual void StartDacq (`
 `int32_t timeout,`
 `int32_t numSubmittedUsbBuffers,`
 `int32_t numUsbBuffers,`
 `int32_t packetsInUrb) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

11.50.3.96 StartDacq() [4/4] `virtual void StartDacq (`
 `int32_t timeout,`
 `int32_t numSubmittedUsbBuffers,`
 `int32_t numUsbBuffers,`
 `int32_t packetsInUrb,`
 `uint32_t virtualDevice) [virtual]`

Start the data acquisition thread and sampling.

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

11.50.3.97 StartLoop() [1/4] `virtual void StartLoop () [virtual]`

Start the data acquisition thread.

11.50.3.98 StartLoop() [2/4] `virtual void StartLoop (
int32_t timeout) [virtual]`

Start the data acquisition thread.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

11.50.3.99 StartLoop() [3/4] `virtual void StartLoop (
int32_t timeout,
int32_t numSubmittedUsbBuffers,
int32_t numUsbBuffers,
int32_t packetsInUrb) [virtual]`

Start the data acquisition thread.

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

11.50.3.100 StartLoop() [4/4] virtual void StartLoop (
 int32_t timeout,
 int32_t numSubmittedUsbBuffers,
 int32_t numUsbBuffers,
 int32_t packetsInUrb,
 uint32_t virtualDevice) [virtual]

Start the data acquisition thread.

Parameters

<i>numSubmittedUsbBuffers</i>	Number of USB Buffers that are simultaneously submitted.
-------------------------------	--

Parameters

<i>timeout</i>	Timeout in ms.
----------------	----------------

Parameters

<i>numUsbBuffers</i>	Number of USB Buffers to use.
----------------------	-------------------------------

Parameters

<i>packetsInUrb</i>	Packets in each URB.
---------------------	----------------------

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

11.50.3.101 StopDacq() [1/2] `virtual void StopDacq () [virtual]`

Stop the data acquisition thread and sampling.

11.50.3.102 StopDacq() [2/2] `virtual void StopDacq (
uint32_t virtualDevice) [virtual]`

Stop the data acquisition thread and sampling.

Parameters

<i>virtualDevice</i>	Virtual Device to start.
----------------------	--------------------------

11.50.3.103 StopLoop() `virtual void StopLoop () [virtual]`

11.50.4 Member Data Documentation

11.50.4.1 Error_Callback_Aquisition_Stopped `const int Error_Callback_Aquisition_Stopped = 0x200 [static]`

11.50.4.2 Error_Callback_Data_lost `const int Error_Callback_Data_lost = 5 [static]`

11.50.4.3 Error_Callback_Frames_Lost `const int Error_Callback_Frames_Lost = 4 [static]`

11.50.4.4 Error_Callback_Packet_Error `const int Error_Callback_Packet_Error = 1 [static]`

11.50.4.5 Error_Callback_Queue_Full `const int Error_Callback_Queue_Full = 0x100 [static]`

11.50.4.6 Error_Callback_RingQueue_Full `const int Error_Callback_RingQueue_Full = 3 [static]`

11.50.5 Property Documentation

11.50.5.1 Samplerate `virtual int Samplerate [get], [set]`

The sampling frequency of the device in Hz.


11.50.6 Event Documentation

11.50.6.1 ChannelDataEvent `OnChannelData^ ChannelDataEvent [add], [remove], [raise]`

11.50.6.2 ErrorEvent `OnError^ ErrorEvent [add], [remove], [raise]`

11.51 CMcsUsbDeviceStatePushFunctionNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushFunctionNet:

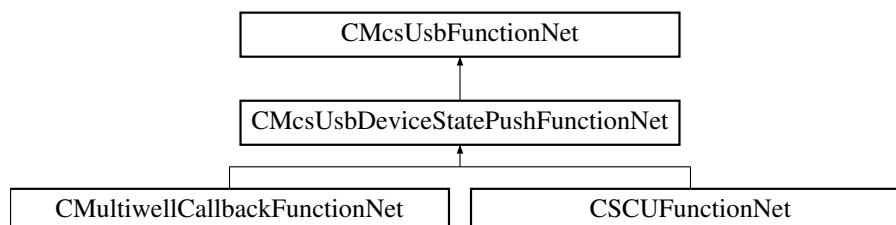

Public Member Functions

- void [TriggerStatus](#) ()

Protected Member Functions

- `CMcsUsbDeviceStatePushFunctionNet` (`CMcsUsbNet`[^] *mcsusb*, `CMcsUsbFunctionPointerContainer`[^] *pDevice*)

Events

- `OnMcsUsbDeviceState`[^] `McsUsbDeviceStateEvent` [add, remove, raise]

Additional Inherited Members

11.51.1 Constructor & Destructor Documentation

11.51.1.1 `CMcsUsbDeviceStatePushFunctionNet()` `CMcsUsbDeviceStatePushFunctionNet` (
`CMcsUsbNet`[^] *mcsusb*,
`CMcsUsbFunctionPointerContainer`[^] *pDevice*) [protected]

11.51.2 Member Function Documentation


11.51.2.1 `TriggerStatus()` `void TriggerStatus` ()

11.51.3 Event Documentation

11.51.3.1 `McsUsbDeviceStateEvent` `OnMcsUsbDeviceState`[^] `McsUsbDeviceStateEvent` [add], [remove], [raise]

11.52 `CMcsUsbDeviceStatePushNet` Class Reference

Inheritance diagram for `CMcsUsbDeviceStatePushNet`:

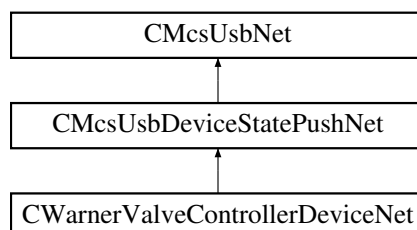

Public Member Functions

- void [TriggerStatus](#) ()

Protected Member Functions

- [CMcsUsbDeviceStatePushNet](#) ([CMcsUsbPointerContainer](#)^ pDevice)

Events

- [OnMcsUsbDeviceState](#)^ [McsUsbDeviceStateEvent](#) [add, remove, raise]

Additional Inherited Members**11.52.1 Constructor & Destructor Documentation**

11.52.1.1 CMcsUsbDeviceStatePushNet() [CMcsUsbDeviceStatePushNet](#) (
[CMcsUsbPointerContainer](#)^ pDevice) [protected]

11.52.2 Member Function Documentation


11.52.2.1 TriggerStatus() void [TriggerStatus](#) ()

11.52.3 Event Documentation

11.52.3.1 McsUsbDeviceStateEvent [OnMcsUsbDeviceState](#)^ [McsUsbDeviceStateEvent](#) [add], [remove], [raise]

11.53 CMcsUsbFactoryNet Class Reference

Inheritance diagram for CMcsUsbFactoryNet:

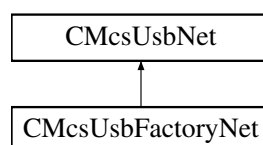

Public Member Functions

- [CMcsUsbFactoryNet](#) ()
- [~CMcsUsbFactoryNet](#) ()
- unsigned int [GetNumDestinations](#) ()
- String ^ [GetDestinationName](#) (unsigned int index)
- String ^ [GetDestinationName](#) (CFirmwareDestinationNet dest)
- void [SetDestinationSerialNumber](#) (CFirmwareDestinationNet dest, String^ serialNumber)
- String ^ [GetDestinationSerialNumber](#) (CFirmwareDestinationNet dest)
- CFirmwareDestinationNet [GetDestination](#) (unsigned int index)
- CFirmwareDestinationNet [GetDestination](#) (String^ Key)
- unsigned int [GetDestinationTargetAddress](#) (CFirmwareDestinationNet destination)
Gets the target base address for the destination.
- uint32_t [ChangeSerialNumber](#) (String^ serial)
- bool [LoadUserFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry)
Send the DSP Firmware to the MEA21 device.
- bool [LoadUserFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, uint32_t LockMask)
- bool [UpdateFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange^ deleg, OnUpdateFirmwareProgress^ progress, bool SkipWait)
Flashes a firmware file to the device.
- bool [UpdateFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange^ deleg, OnUpdateFirmwareProgress^ progress, bool SkipWait, unsigned int LockMask)
- bool [UpdateFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet dest)
Flashes a firmware file to the device.
- bool [UpdateFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet dest, bool SkipWait)
Flashes a firmware file to the device.
- bool [UpdateFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet dest, bool SkipWait, uint32_t LockMask)
- bool [CompareFirmware](#) (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange^ deleg, OnUpdateFirmwareProgress^ progress, String^ MessagePrefix, unsigned int LockMask, [System::Runtime::InteropServices::Out] String^% ErrorText, [System::Runtime::InteropServices::Out] String^% Protokoll)
- uint32_t [Coldstart](#) (CFirmwareDestinationNet dest)
- int32_t [GetXilinxFlashOffset](#) (CFirmwareDestinationNet dest)
- uint32_t [GetXilinxFlashReadCommand](#) (CFirmwareDestinationNet dest)
- array< uint8_t > ^ [DownloadFirmware](#) (CFirmwareDestinationNet Dest, uint32_t Address, uint32_t length)
- bool [GetUserCodeFromFlash](#) (unsigned int FPGA, unsigned int Address, [System::Runtime::InteropServices::Out] unsigned int% Usercode)
- array< unsigned char > ^ [ReadBlockFromFlash](#) (unsigned int FPGA, unsigned int Address)
- void [ReadBlockFromFlash](#) (unsigned int FPGA, unsigned int Address, array< unsigned char >^ buffer, int position)
- array< unsigned char > ^ [ReadBlockFromIFBGlobalEEPROM](#) (unsigned int Address)
- array< unsigned char > ^ [ReadBlockFromNVMEM](#) (unsigned int FPGA, unsigned int Offset, unsigned int Address)

Static Public Member Functions

- static String ^ [GetDestinationDisplayLabel](#) (String^ RawLabel, [CFirmwareDestinationNet](#) dest)
- static String ^ [FindFirmwareVersionMagicInBuffer](#) (array< unsigned char >^ buffer, int length, [System::Runtime::InteropServices::Out] int% position)
- static bool [GetFirmwareVersionFromFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
Retrives version info from a Firmware update file.
- static bool [GetFirmwareVersionFromFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version, [System::Runtime::InteropServices::Out] uint32_t% Position)
- static bool [GetFirmwareVersionFromHexFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
- static uint32_t [GetChecksumFromFX3Image](#) (String^ FirmwareFile)
- static uint32_t [GetUSBDeviceIDFromFX3Image](#) (String^ FirmwareFile)
- static bool [GetUserCodeFromBitFile](#) (String^ FirmwareFile, [System::Runtime::InteropServices::Out] unsigned int% Usercode)

Static Public Attributes

- static const uint32_t [FX3MCSDDataAddress](#) = 0x40037E00
- static const uint32_t [FX3MCSDDataDeviceIdOffset](#) = 0x4
- static const uint32_t [FX3MCSDDataVersionOffset](#) = 0x8
- static const uint32_t [FX3MCSDDataIFB2ImageOffset](#) = 0xC
- static const uint32_t [FX3MCSDDataIFB1ImageOffset](#) = 0x2C

Additional Inherited Members

11.53.1 Constructor & Destructor Documentation

11.53.1.1 [CMcsUsbFactoryNet\(\)](#) [CMcsUsbFactoryNet](#) ()

11.53.1.2 [~CMcsUsbFactoryNet\(\)](#) [~CMcsUsbFactoryNet](#) ()

11.53.2 Member Function Documentation

11.53.2.1 [ChangeSerialNumber\(\)](#) uint32_t ChangeSerialNumber (String^ serial)

11.53.2.2 Coldstart() uint32_t Coldstart (
 CFirmwareDestinationNet dest)

11.53.2.3 CompareFirmware() bool CompareFirmware (
 String^ FirmwareFile,
 CMcsUsbListEntryNet^ listEntry,
 CFirmwareDestinationNet Dest,
 OnUpdateFirmwareStatusChange^ deleg,
 OnUpdateFirmwareProgress^ progress,
 String^ MessagePrefix,
 unsigned int LockMask,
 [System::Runtime::InteropServices::Out] String^% ErrorText,
 [System::Runtime::InteropServices::Out] String^% Protokoll)

11.53.2.4 DownloadFirmware() array<uint8_t> ^ DownloadFirmware (
 CFirmwareDestinationNet Dest,
 uint32_t Address,
 uint32_t length)

11.53.2.5 FindFirmwareVersionMagicInBuffer() static String ^ FindFirmwareVersionMagicInBuffer (
 array< unsigned char >^ buffer,
 int length,
 [System::Runtime::InteropServices::Out] int% position) [static]

11.53.2.6 GetChecksumFromFX3Image() static uint32_t GetChecksumFromFX3Image (
 String^ FirmwareFile) [static]

11.53.2.7 GetDestination() [1/2] CFirmwareDestinationNet GetDestination (
 String^ Key)

11.53.2.8 GetDestination() [2/2] CFirmwareDestinationNet GetDestination (
 unsigned int index)

11.53.2.9 GetDestinationDisplayLabel() static String ^ GetDestinationDisplayLabel (
 String^ RawLabel,
 CFirmwareDestinationNet dest) [static]

11.53.2.10 GetDestinationName() [1/2] `String ^ GetDestinationName (CFirmwareDestinationNet dest)`

11.53.2.11 GetDestinationName() [2/2] `String ^ GetDestinationName (unsigned int index)`

11.53.2.12 GetDestinationSerialNumber() `String ^ GetDestinationSerialNumber (CFirmwareDestinationNet dest)`

11.53.2.13 GetDestinationTargetAddress() `unsigned int GetDestinationTargetAddress (CFirmwareDestinationNet destination)`

Gets the target base address for the destination.

Parameters

<i>destination</i>	The destination to be queried.
--------------------	--------------------------------

Returns

The base address as a 32 bit number, only the lower 16 bit represent the address.

11.53.2.14 GetFirmwareVersionFromFile() [1/2] `static bool GetFirmwareVersionFromFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version) [static]`

Retrives version info from a Firmware update file.

11.53.2.15 GetFirmwareVersionFromFile() [2/2] `static bool GetFirmwareVersionFromFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version, [System::Runtime::InteropServices::Out] uint32_t% Position) [static]`

11.53.2.16 GetFirmwareVersionFromHexFile() `static bool GetFirmwareVersionFromHexFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version) [static]`

11.53.2.17 GetNumDestinations() `unsigned int GetNumDestinations ()`

11.53.2.18 GetUSBDeviceIDFromFX3Image() `static uint32_t GetUSBDeviceIDFromFX3Image (String^ FirmwareFile) [static]`

11.53.2.19 GetUsercodeFromBitFile() `static bool GetUsercodeFromBitFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] unsigned int% Usercode) [static]`

11.53.2.20 GetUsercodeFromFlash() `bool GetUsercodeFromFlash (unsigned int FPGA, unsigned int Address, [System::Runtime::InteropServices::Out] unsigned int% Usercode)`

11.53.2.21 GetXilinxFlashOffset() `int32_t GetXilinxFlashOffset (CFirmwareDestinationNet dest)`

11.53.2.22 GetXilinxFlashReadCommand() `uint32_t GetXilinxFlashReadCommand (CFirmwareDestinationNet dest)`

11.53.2.23 LoadUserFirmware() `[1/2] bool LoadUserFirmware (String^ FirmwareFile, CMcsUsbListEntryNet^ listEntry)`

Send the DSP Firmware to the MEA21 device.

Parameters

<i>FirmwareFile</i>	Filename of the DSP Firmware (*.bin) file.
---------------------	--

Parameters

<i>listEntry</i>	Device to use for the connection. See CMcsUsbListNet .
------------------	--

11.53.2.24 LoadUserFirmware() [2/2] `bool LoadUserFirmware (`
 `String^ FirmwareFile,`
 `CMcsUsbListEntryNet^ listEntry,`
 `uint32_t LockMask)`

11.53.2.25 ReadBlockFromFlash() [1/2] `array<unsigned char> ^ ReadBlockFromFlash (`
 `unsigned int FPGA,`
 `unsigned int Address)`

11.53.2.26 ReadBlockFromFlash() [2/2] `void ReadBlockFromFlash (`
 `unsigned int FPGA,`
 `unsigned int Address,`
 `array< unsigned char >^ buffer,`
 `int position)`

11.53.2.27 ReadBlockFromIFBGlobalEEProm() `array<unsigned char> ^ ReadBlockFromIFBGlobal↵`
`EEProm (`
 `unsigned int Address)`

11.53.2.28 ReadBlockFromNVMEM() `array<unsigned char> ^ ReadBlockFromNVMEM (`
 `unsigned int FPGA,`
 `unsigned int Offset,`
 `unsigned int Address)`

11.53.2.29 SetDestinationSerialNumber() `void SetDestinationSerialNumber (`
 `CFirmwareDestinationNet dest,`
 `String^ serialnumber)`

11.53.2.30 UpdateFirmware() [1/5] `bool UpdateFirmware (`
 `String^ FirmwareFile,`
 `CMcsUsbListEntryNet^ listEntry,`
 `CFirmwareDestinationNet dest)`

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------

Parameters

<i>listEntry</i>	Device to use for the connection.
------------------	-----------------------------------

11.53.2.31 UpdateFirmware() [2/5] `bool UpdateFirmware (`
 `String^ FirmwareFile,`
 `CMcsUsbListEntryNet^ listEntry,`
 `CFirmwareDestinationNet dest,`
 `bool SkipWait)`

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------

Parameters

<i>listEntry</i>	Device to use for the connection.
------------------	-----------------------------------

11.53.2.32 UpdateFirmware() [3/5] `bool UpdateFirmware (`
 `String^ FirmwareFile,`
 `CMcsUsbListEntryNet^ listEntry,`
 `CFirmwareDestinationNet dest,`
 `bool SkipWait,`
 `uint32_t LockMask)`

11.53.2.33 UpdateFirmware() [4/5] `bool UpdateFirmware (`
 `String^ FirmwareFile,`
 `CMcsUsbListEntryNet^ listEntry,`
 `CFirmwareDestinationNet Dest,`

```

    OnUpdateFirmwareStatusChange^ deleg,
    OnUpdateFirmwareProgress^ progress,
    bool SkipWait )

```

Flashes a firmware file to the device.

Parameters

<i>FirmwareFile</i>	Filename of the Firmware file.
---------------------	--------------------------------

```

11.53.2.34 UpdateFirmware() [5/5] bool UpdateFirmware (
    String^ FirmwareFile,
    CMcsUsbListEntryNet^ listEntry,
    CFirmwareDestinationNet Dest,
    OnUpdateFirmwareStatusChange^ deleg,
    OnUpdateFirmwareProgress^ progress,
    bool SkipWait,
    unsigned int LockMask )

```

11.53.3 Member Data Documentation

11.53.3.1 FX3MCSDDataAddress const uint32_t FX3MCSDDataAddress = 0x40037E00 [static]

11.53.3.2 FX3MCSDDataDeviceIdOffset const uint32_t FX3MCSDDataDeviceIdOffset = 0x4 [static]

11.53.3.3 FX3MCSDDataIFB1ImageOffset const uint32_t FX3MCSDDataIFB1ImageOffset = 0x2C [static]

11.53.3.4 FX3MCSDDataIFB2ImageOffset const uint32_t FX3MCSDDataIFB2ImageOffset = 0xC [static]

11.53.3.5 FX3MCSDDataVersionOffset const uint32_t FX3MCSDDataVersionOffset = 0x8 [static]

[illegible]

- **CMcsUsbFunctionNet** (**CMcsUsbNet**[^] mcsusb)
- virtual ~**CMcsUsbFunctionNet** (void)
- **ICMcsUsbFunctionNet** ()
- void **ThrowCUsbExceptionNetOnError** (uint32 t status)

- `CMcsUsbFunctionNet` (`CMcsUsbNet`^{mcsusb}, `CMcsUsbFunctionPointerContainer`^{mcsusbfunction})

- CMcsUsbNet ^ m_pMcsUsb
- CMcsUsbFunction * m_pMcsUsbFunction

Generated by Doxygen

11.54.1.1 CMcsUsbFunctionNet() [1/2] `CMcsUsbFunctionNet (`
`CMcsUsbNet^ mcsusb)`

11.54.1.2 ~CMcsUsbFunctionNet() `virtual ~CMcsUsbFunctionNet (`
`void) [virtual]`

11.54.1.3 "!CMcsUsbFunctionNet() `!CMcsUsbFunctionNet ()`

11.54.1.4 CMcsUsbFunctionNet() [2/2] `CMcsUsbFunctionNet (`
`CMcsUsbNet^ mcsusb,`
`CMcsUsbFunctionPointerContainer^ mcsusbfunction) [protected]`

11.54.2 Member Function Documentation

11.54.2.1 ThrowCUsbExceptionNetOnError() `void ThrowCUsbExceptionNetOnError (`
`uint32_t status)`

11.54.3 Member Data Documentation

11.54.3.1 m_pMcsUsb `CMcsUsbNet ^ m_pMcsUsb [protected]`

11.54.3.2 m_pMcsUsbFunction `CMcsUsbFunction* m_pMcsUsbFunction [protected]`

11.55 CMcsUsbFunctionPointerContainer Class Reference

11.56 CMcsUsbListEntryNet Class Reference

McsUsbListEntryNet identifies a connected device.

Public Member Functions

- `~CMcsUsbListEntryNet ()`
- virtual bool `Equals (Object^ obj)` override
Checks weather two `CMcsUsbListEntryNet` represent the same USB device.
- void `SetStringFormat (String ^ format)`
Specify the text the `CMcsUsbListEntryNet.ToString()` function should return. The special code `N` expands to the device name and `S` expands to the serial number of the device.
- virtual String ^ `ToString ()` override

Static Public Member Functions

- static `CMcsUsbListEntryNet ^ GetEntry ()`
Returns one `CMcsUsbListEntryNet` from the list of USB Devices connected to the computer.
- static `CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice)`
Returns one `CMcsUsbListEntryNet` from the list of USB Devices connected to the computer.
- static `CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice, unsigned int index)`
Returns one `CMcsUsbListEntryNet` from the list of USB Devices connected to the computer.
- static unsigned int `GetEntryCount ()`
Returns the number of devices connected to the computer.
- static unsigned int `GetEntryCount (DeviceEnumNet McsUsbDevice)`
Returns the number of devices connected to the computer.

Properties

- String^ `Manufacturer` [get]
The Manufacturer ID of the device represented by this `CMcsUsbListEntryNet`.
- String^ `Product` [get]
The Product ID of the device represented by this `CMcsUsbListEntryNet`.
- String^ `DeviceName` [get]
The device name of the device represented by this `CMcsUsbListEntryNet`.
- String^ `SerialNumber` [get]
The serial number of the device represented by this `CMcsUsbListEntryNet`.
- String^ `HwVersion` [get]
The hardware revision of the device represented by this `CMcsUsbListEntryNet`.
- `DeviceldNet^ Deviceld` [get]

11.56.1 Detailed Description

McsUsbListEntryNet identifies a connected device.

11.56.2 Constructor & Destructor Documentation

11.56.2.1 `~CMcsUsbListEntryNet()` `~CMcsUsbListEntryNet ()`

11.56.3 Member Function Documentation

11.56.3.1 Equals() `virtual bool Equals (
 Object^ obj) [override], [virtual]`

Checks weather two [CMcsUsbListEntryNet](#) represent the same USB device.

Parameters

<i>obj</i>	The CMcsUsbListEntryNet to compare with.
------------	--

11.56.3.2 GetEntry() [1/3] `static CMcsUsbListEntryNet ^ GetEntry () [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

11.56.3.3 GetEntry() [2/3] `static CMcsUsbListEntryNet ^ GetEntry (
 DeviceEnumNet McsUsbDevice) [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

11.56.3.4 GetEntry() [3/3] `static CMcsUsbListEntryNet ^ GetEntry (
 DeviceEnumNet McsUsbDevice,
 unsigned int index) [static]`

Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Parameters

<i>index</i>	number of the entry to use.
--------------	-----------------------------

Returns

A [CMcsUsbListEntryNet](#) to be used to connect to the device.

11.56.3.5 GetEntryCount() [1/2] `static unsigned int GetEntryCount () [static]`

Returns the number of devices connected to the computer.

Returns

The number of devices.

11.56.3.6 GetEntryCount() [2/2] `static unsigned int GetEntryCount (
DeviceEnumNet McsUsbDevice) [static]`

Returns the number of devices connected to the computer.

Parameters

<i>McsUsbDevice</i>	Specifies the type of devices to look for.
---------------------	--

Returns

The number of devices.

11.56.3.7 SetStringFormat() `void SetStringFormat (
String ^ format)`

Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

<i>format</i>	A String containing the format template. Default is "%N (%S)".
---------------	--

11.56.3.8 ToString() `virtual String ^ ToString () [override], [virtual]`

11.56.4 Property Documentation

11.56.4.1 DeviceId `DeviceIdNet^ DeviceId [get]`

11.56.4.2 DeviceName `String^ DeviceName [get]`

The device name of the device represented by this [CMcsUsbListEntryNet](#).

11.56.4.3 HwVersion `String^ HwVersion [get]`

The hardware revision of the device represented by this [CMcsUsbListEntryNet](#).

11.56.4.4 Manufacturer `String^ Manufacturer [get]`

The Manufacturer ID of the device represented by this [CMcsUsbListEntryNet](#).

11.56.4.5 Product `String^ Product [get]`

The Product ID of the device represented by this [CMcsUsbListEntryNet](#).

11.56.4.6 SerialNumber `String^ SerialNumber [get]`

The serial number of the device represented by this [CMcsUsbListEntryNet](#).

11.57 CMcsUsbListNet Class Reference

Class to handle a list of connected MCS USB devices.

Public Member Functions

- [CMcsUsbListNet](#) ([DeviceEnumNet](#) McsUsbDevice)
Initializes a new instance of [CMcsUsbListNet](#) class.
- [CMcsUsbListNet](#) (array< [DeviceIdNet](#)[^]> ^ DeviceIdList)
Initializes a new instance of [CMcsUsbListNet](#) class.
- [~CMcsUsbListNet](#) ()
Destructor: called by [Dispose\(\)](#)
- [!CMcsUsbListNet](#) ()
Finalizer: called by GC before collecting
- void [SetStringFormat](#) (String ^ format)
Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.
- [uint32_t](#) [GetNumberOfDevices](#) ()
Gets the number of devices currently in the list.
- [CMcsUsbListEntryNet](#)[^] [GetUsbListEntry](#) (unsigned int index)
Returns one [CMcsUsbListEntryNet](#) from the list of USB Devices connected to the computer.
- array< [CMcsUsbListEntryNet](#)[^]> ^ [GetUsbListEntries](#) ()
Returns all entries from the list of USB Devices connected to the computer.
- bool [IsDeviceTypeOf](#) ([CMcsUsbListEntryNet](#)[^] entry, [DeviceEnumNet](#) McsUsbDevice)

Properties

- [uint32_t](#) [Count](#) [get]
Gets the number of devices currently in the list.

Events

- [OnDeviceArrivalRemoval](#)[^] [DeviceArrival](#)
- [OnDeviceArrivalRemoval](#)[^] [DeviceRemoval](#)

11.57.1 Detailed Description

Class to handle a list of connected MCS USB devices.

11.57.2 Constructor & Destructor Documentation

11.57.2.1 CMcsUsbListNet() [1/2] [CMcsUsbListNet](#) ([DeviceEnumNet](#) McsUsbDevice)

Initializes a new instance of [CMcsUsbListNet](#) class.

11.57.2.2 CMcsUsbListNet() [2/2] `CMcsUsbListNet (array< DeviceIdNet^>^ DeviceIdList)`

Initializes a new instance of `CMcsUsbListNet` class.

11.57.2.3 ~CMcsUsbListNet() `~CMcsUsbListNet ()`

Destructor: called by `Dispose()`

11.57.2.4 ~!CMcsUsbListNet() `!CMcsUsbListNet ()`

Finalizer: called by GC before collecting

11.57.3 Member Function Documentation

11.57.3.1 GetNumberOfDevices() `uint32_t GetNumberOfDevices ()`

Gets the number of devices currently in the list.

Returns

The number of devices currently in the list.

11.57.3.2 GetUsbListEntries() `array<CMcsUsbListEntryNet^> ^ GetUsbListEntries ()`

Returns all entries from the list of USB Devices connected to the computer.

11.57.3.3 GetUsbListEntry() `CMcsUsbListEntryNet ^ GetUsbListEntry (unsigned int index)`

Returns one `CMcsUsbListEntryNet` from the list of USB Devices connected to the computer.

Parameters

<i>index</i>	number of the entry to use.
--------------	-----------------------------

11.57.3.4 IsDeviceTypeOf() `bool IsDeviceTypeOf (`
 `CMcsUsbListEntryNet^ entry,`
 `DeviceEnumNet McsUsbDevice)`

11.57.3.5 SetStringFormat() `void SetStringFormat (`
 `String ^ format)`

Specify the text the [CMcsUsbListEntryNet.ToString\(\)](#) function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

<i>format</i>	A String containing the format template. Default is "%N (%S)".
---------------	--

11.57.4 Property Documentation

11.57.4.1 Count `uint32_t Count [get]`

Gets the number of devices currently in the list.

11.57.5 Event Documentation


11.57.5.1 DeviceArrival `OnDeviceArrivalRemoval^ DeviceArrival`

11.57.5.2 DeviceRemoval `OnDeviceArrivalRemoval^ DeviceRemoval`

11.58 CMcsUsbNet Class Reference

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

Inheritance diagram for CMcsUsbNet:

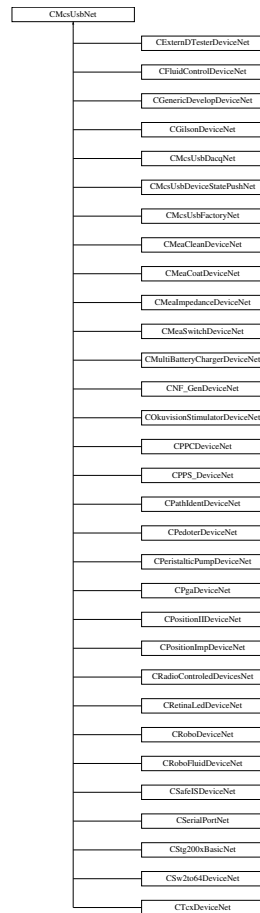

Public Member Functions

- [CMcsUsbNet](#) ()
Initializes a new instance of the base class to handle MCS USB devices.
- [CMcsUsbNet](#) ([McsBusTypeEnumNet](#) bustype)
Initializes a new instance of the base class to handle MCS USB devices.
- virtual [~CMcsUsbNet](#) ()
- [!CMcsUsbNet](#) ()
- [DeviceEnumNet](#) [GetDeviceEnum](#) ()
- virtual uint32_t [Connect](#) ()
Opens a connection to the device.
- virtual uint32_t [Connect](#) (unsigned int LockMask)
Opens a connection to the device.
- virtual uint32_t [Connect](#) ([CMcsUsbListEntryNet](#)[^] entry)
Opens a connection to the device.
- virtual uint32_t [Connect](#) ([CMcsUsbListEntryNet](#)[^] entry, unsigned int LockMask)
Opens a connection to the device.

- virtual uint32_t [GetStatus](#) ([System::Runtime::InteropServices::Out]uint32_t% iStatus)
- virtual bool [IsConnected](#) ()
Check if a device is Connected.
- virtual void [Disconnect](#) ()
Disconnect from a device.
- [CMcsUsbListEntryNet](#) ^ [GetUsbListEntry](#) ()
- virtual String ^ [GetSerialNumber](#) ()
Query the Serial Number of the device.
- [DriverVersionNet](#) ^ [GetVersion](#) ()
- [DriverVersionNet](#) ^ [GetVersion](#) (CFirmwareDestinationNet dest)
- [DeviceldNet](#) ^ [GetDeviceld](#) ()
- void [MultibootSelectImage](#) (unsigned int sector)
Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.
- String ^ [MultibootGetImageld](#) (unsigned int sector)
Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.
- uint32_t [MultibootGetCypressImageld](#) (unsigned int sector)
Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.
- uint32_t [MultibootGetSelectedImage](#) ()
Gets sector index of selected FPGA boot image on IFB
- uint32_t [GetMea21UsbPort](#) ()
Gets the USB port if an IFB that is used by this connection
- [HeadstageIdEnumNet](#) [GetHeadstageID](#) (uint32_t headstage)
Gets the ID of a connected headstage.
- bool [GetHeadstagePresent](#) (uint32_t headstage)
queries whether a headstage is present
- bool [GetHeadstageActive](#) (uint32_t headstage)
queries whether a headstage is active
- void [RescanHeadstage](#) (uint32_t headstage)
rescans and activates a headstage
- array< BYTE > ^ [GetSoftwareKey](#) (unsigned int index)
- void [SetSoftwareKey](#) (unsigned int index, array< BYTE >^ buffer)
- void [RemoveSoftwareKey](#) (unsigned int index)
- void [AddSoftwareKey](#) (String^ key)
- bool [EmptyKey](#) (String^ key)
- bool [ValidKey](#) (String^ key, [System::Runtime::InteropServices::Out]String^% serial_number)
- bool [ValidKey](#) (String^ key, uint8_t ProgrammID, uint8_t majorversion, [System::Runtime::InteropServices::Out]String^% serial_number)
- bool [HasSoftwareKey](#) (uint8_t ProgrammID, uint8_t majorversion)
- bool [HasSoftwareKey](#) (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- String ^ [GetSoftwareKeyString](#) (uint8_t ProgrammID, uint8_t majorversion)
- String ^ [GetSoftwareKeyString](#) (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- bool [IsDeviceHighSpeedCapable](#) ()
- bool [IsDeviceHighSpeed](#) ()
- [McsUsbSpeedEnumNet](#) [GetDeviceCapableSpeed](#) ()
- [McsUsbSpeedEnumNet](#) [GetDeviceSpeed](#) ()
Query the Connection Speed of the device.
- unsigned int [TxnTestMemoryWrite](#) (unsigned short index)
- unsigned int [TxnTestMemoryReadAndCheck](#) (unsigned short index)
- void [TxnSetSerialNumber](#) (unsigned int number)
- unsigned int [TxnGetSerialNumber](#) ()
- unsigned int [ReadRegister](#) (unsigned int reg)
- array< uint32_t > ^ [ReadRegister](#) (unsigned int reg, int length)

- unsigned int [ReadRegister32](#) (unsigned int adr)
- unsigned int [ReadRegisterTimeSlot](#) (unsigned int reg, int TimeSlot)
- void [WriteRegister](#) (unsigned int reg, unsigned int value)
- void [WriteRegisterValue](#) (unsigned int reg, unsigned int value)
- void [WriteRegister32](#) (unsigned int adr, unsigned int value)
- void [WriteRegister](#) (unsigned int reg, array< unsigned int >^ values)
- void [WriteRegisterArray](#) (unsigned int reg, array< unsigned int >^ values)
- void [WriteRegisterTimeSlot](#) (unsigned int reg, unsigned int value, int TimeSlot)
- void [WriteRegisterTimeSlot](#) (unsigned int reg, array< unsigned int >^ values, int TimeSlot)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t DMA_value)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t DMA_value, uint32_t EEPROMSize)
- bool [ReadEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::InteropServices::Out]uint32_t DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize)
- void [WriteEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize)
- void [EraseEepromRegisterPreconfig](#) (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- unsigned int [GetLastUSBError](#) ()
- void [ThrowCUsbExceptionNetOnError](#) (uint32_t status)
- bool [GetDeviceCannotStallOutRequests](#) ()
- String ^ [GetHardwareRevision](#) ()
- unsigned int [GetFirmwareVersion](#) (CFirmwareDestinationNet destination)
Gets the firmware version for the destination.
- uint8_t [GetNumConfigurations](#) ()
- uint8_t [GetConfiguration](#) ()
- void [SetConfiguration](#) (uint8_t config)
- uint32_t [GetDeviceRootHubVendorID](#) ()
Gets the Vendor ID of the USB root hub the device is connected to.
- [UsbVendorIdEnumNet GetDeviceRootHubVendorEnum](#) ()
Gets the Vendor ID of the USB root hub the device is connected to.
- String ^ [GetDeviceRootHubVendorName](#) ()
Gets the Vendor Name of the USB root hub the device is connected to.
- void [EnableExceptions](#) (bool enable)
Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.
- bool [IsExceptionsEnabled](#) ()
- uint32_t [GetStatusOfLastCommand](#) ()
Gets the status of the last call to the McsUsb Library.
- uint32_t [CyclePort](#) ()
- void [AssociateToThis](#) (CMcsUsbNet^ device)

Static Public Member Functions

- static String ^ [GetErrorText](#) (unsigned int Status)
Gets the error text string that belongs to a status number.

Static Public Attributes

- static const uint32_t [Status_Crc](#) = (0xE0100001L)
- static const uint32_t [Status_Btstuff](#) = (0xE0100002L)
- static const uint32_t [Status_DataToggleMismatch](#) = (0xE0100003L)
- static const uint32_t [Status_Stall](#) = (0xE0100004L)
- static const uint32_t [Status_DevNotResponding](#) = (0xE0100005L)
- static const uint32_t [Status_PidCheckFailure](#) = (0xE0100006L)
- static const uint32_t [Status_UnexpectedPid](#) = (0xE0100007L)
- static const uint32_t [Status_DataOverrun](#) = (0xE0100008L)
- static const uint32_t [Status_DataUnderrun](#) = (0xE0100009L)
- static const uint32_t [Status_BufferOverrun](#) = (0xE010000CL)
- static const uint32_t [Status_BufferUnderrun](#) = (0xE010000DL)
- static const uint32_t [Status_NotAccessed](#) = (0xE010000FL)
- static const uint32_t [Status_Fifo](#) = (0xE0100010L)
- static const uint32_t [Status_EndpointHalted](#) = (0xE0100030L)
- static const uint32_t [Status_NoMemory](#) = (0xE0100100L)
- static const uint32_t [Status_InvalidUrbFunction](#) = (0xE0100200L)
- static const uint32_t [Status_InvalidParameter](#) = (0xE0100300L)
- static const uint32_t [Status_InvalidDeviceHandle](#) = (0xE0100013L)
- static const uint32_t [Status_InvalidHandle](#) = (0xE0100012L)
- static const uint32_t [Status_ErrorBusy](#) = (0xE0100400L)
- static const uint32_t [Status_RequestFailed](#) = (0xE0100500L)
- static const uint32_t [Status_InvalidPipeHandle](#) = (0xE0100600L)
- static const uint32_t [Status_NoBandwidth](#) = (0xE0100700L)
- static const uint32_t [Status_InternalHcError](#) = (0xE0100800L)
- static const uint32_t [Status_ErrorShortTransfer](#) = (0xE0100900L)
- static const uint32_t [Status_BadStartFrame](#) = (0xE0100A00L)
- static const uint32_t [Status_IsochRequestFailed](#) = (0xE0100B00L)
- static const uint32_t [Status_FrameControlOwned](#) = (0xE0100C00L)
- static const uint32_t [Status_ControlNotOwned](#) = (0xE0100D00L)
- static const uint32_t [Status_Canceled](#) = (0xE0110000L)
- static const uint32_t [Status_Canceling](#) = (0xE0120000L)
- static const uint32_t [Status_AlreadyConfigured](#) = (0xE0110001L)
- static const uint32_t [Status_Unconfigured](#) = (0xE0110002L)
- static const uint32_t [Status_NoSuchDevice](#) = (0xE01F0002L)
- static const uint32_t [Status_DeviceNotFound](#) = (0xE01F0003L)
- static const uint32_t [Status_NotSupported](#) = (0xE01F0005L)
- static const uint32_t [Status_IoPending](#) = (0xE01F0006L)
- static const uint32_t [Status_IoTimeout](#) = (0xE01F0007L)
- static const uint32_t [Status_DeviceRemoved](#) = (0xE01F0008L)
- static const uint32_t [Status_PipeNotLinked](#) = (0xE01F0009L)
- static const uint32_t [Status_ConnectedPipes](#) = (0xE01F000AL)
- static const uint32_t [Status_DeviceLocked](#) = (0xE01F0010L)
- static const uint32_t [Status_RequestMutexTimeout](#) = (0xE01F0020L)
- static const uint32_t [Status_RequestMutexFailed](#) = (0xE01F0021L)
- static const uint32_t [Status_LastUsbErrorMismatch](#) = (0xE01F0022L)
- static const uint32_t [WPAError_ScanningsPending](#) = ((0xA0220000L) | 0x0036)

Properties

- virtual String^ [SerialNumber](#) [get]

11.58.1 Detailed Description

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

11.58.2 Constructor & Destructor Documentation

11.58.2.1 CMcsUsbNet() [1/2] `CMcsUsbNet ()`

Initializes a new instance of the base class to handle MCS USB devices.

11.58.2.2 CMcsUsbNet() [2/2] `CMcsUsbNet (McsBusTypeEnumNet bustype)`

Initializes a new instance of the base class to handle MCS USB devices.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

11.58.2.3 ~CMcsUsbNet() `virtual ~CMcsUsbNet () [virtual]`

11.58.2.4 "!CMcsUsbNet() `!CMcsUsbNet ()`

11.58.3 Member Function Documentation

11.58.3.1 AddSoftwareKey() `void AddSoftwareKey (String^ key)`

11.58.3.2 AssociateToThis() `void AssociateToThis (CMcsUsbNet^ device)`

11.58.3.3 Connect() [1/4] virtual uint32_t Connect () [virtual]

Opens a connection to the device.

Returns

Error Status. 0 on success.

11.58.3.4 Connect() [2/4] virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry) [virtual]

Opens a connection to the device.

Parameters

<i>entry</i>	The Device List Entry for the device to be connected.
--------------	---

Returns

Error Status. 0 on success.

11.58.3.5 Connect() [3/4] virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry,
unsigned int LockMask) [virtual]

Opens a connection to the device.

Parameters

<i>entry</i>	The Device List Entry for the device to be connected.
<i>LockMask</i>	The Lock Mask for this connection.

Returns

Error Status. 0 on success.

11.58.3.6 Connect() [4/4] virtual uint32_t Connect (
unsigned int LockMask) [virtual]

Opens a connection to the device.

Parameters

<i>LockMask</i>	The Lock Mask for this connection.
-----------------	------------------------------------

Returns

Error Status. 0 on success.

11.58.3.7 CyclePort() `uint32_t CyclePort ()`

11.58.3.8 Disconnect() `virtual void Disconnect () [virtual]`

Disconnect from a device.

11.58.3.9 EmptyKey() `bool EmptyKey (`
 `String^ key)`

11.58.3.10 EnableExceptions() `void EnableExceptions (`
 `bool enable)`

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the `GetStatusOfLastCommand` call instead.

Parameters

<i>enable</i>	True to enable Exceptions, False to disable.
---------------	--

11.58.3.11 EraseEepromRegisterPreconfig() [1/3] `void EraseEepromRegisterPreconfig (`
 `uint32_t EEPROMBase,`
 `uint32_t DMA_reg)`

11.58.3.12 EraseEepromRegisterPreconfig() [2/3] `void EraseEepromRegisterPreconfig (`
 `uint32_t EEPROMBase,`
 `uint32_t DMA_reg,`
 `uint32_t EEPROMSize)`

11.58.3.13 EraseEepromRegisterPreconfig() [3/3] void EraseEepromRegisterPreconfig (
 uint32_t *EEPROMBase*,
 uint32_t *DMA_reg*,
 uint32_t *EEPROMSize*,
 uint32_t *EepromStartAddress*)

11.58.3.14 GetConfiguration() uint8_t GetConfiguration ()

11.58.3.15 GetDeviceCannotStallOutRequests() bool GetDeviceCannotStallOutRequests ()

11.58.3.16 GetDeviceCapableSpeed() McsUsbSpeedEnumNet GetDeviceCapableSpeed ()

11.58.3.17 GetDeviceEnum() DeviceEnumNet GetDeviceEnum ()

11.58.3.18 GetDeviceId() DeviceIdNet ^ GetDeviceId ()

11.58.3.19 GetDeviceRootHubVendorEnum() UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ()

Gets the Vendor ID of the USB root hub the device is connected to.

Returns

An enum which enumerates the PCI Vendor ID.

11.58.3.20 GetDeviceRootHubVendorID() uint32_t GetDeviceRootHubVendorID ()

Gets the Vendor ID of the USB root hub the device is connected to.

Returns

The PCI Vendor ID, 0x8086 for Intel, 0x1912 for Renesas, 0x1b21 for ASMedia.

11.58.3.21 GetDeviceRootHubVendorName() `String ^ GetDeviceRootHubVendorName ()`

Gets the Vendor Name of the USB root hub the device is connected to.

Returns

The PCI Vendor Name, either "Intel", "Renesas", "ASMedia" or "unknown".

11.58.3.22 GetDeviceSpeed() `McsUsbSpeedEnumNet GetDeviceSpeed ()`

Query the Connection Speed of the device.

Returns

0 for Low-Speed, 1 for Full-Speed, 2 for High-Speed and 3 for SuperSpeed.

11.58.3.23 GetErrorText() `static String ^ GetErrorText (unsigned int Status) [static]`

Gets the error text string that belongs to a status number.

Parameters

<i>Status</i>	The status number you want the text for.
---------------	--

Returns

The error text string that belongs to the status number.

11.58.3.24 GetFirmwareVersion() `unsigned int GetFirmwareVersion (CFirmwareDestinationNet destination)`

Gets the firmware version for the destination.

Parameters

<i>destination</i>	The destination to be queried.
--------------------	--------------------------------

Returns

The firmware version as a 32 bit number, the upper 16 bit contain the major version number, the lower 16 bit the minor version number.

11.58.3.25 GetHardwareRevision() `String ^ GetHardwareRevision ()`

11.58.3.26 GetHeadstageActive() `bool GetHeadstageActive (
 uint32_t headstage)`

queries whether a headstage is active

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

true if the headstage is active

11.58.3.27 GetHeadstageID() `HeadstageIdEnumNet GetHeadstageID (
 uint32_t headstage)`

Gets the ID of a connected headstage.

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

enumerated Headstage ID

11.58.3.28 GetHeadstagePresent() `bool GetHeadstagePresent (
 uint32_t headstage)`

queries whether a headstage is present

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

Returns

true if the headstage is present

11.58.3.29 GetLastUSBError() `unsigned int GetLastUSBError ()`

11.58.3.30 GetMea21UsbPort() `uint32_t GetMea21UsbPort ()`

Gets the USB port if an IFB that is used by this connection

Returns

number of used port; range: 0..1

11.58.3.31 GetNumConfigurations() `uint8_t GetNumConfigurations ()`

11.58.3.32 GetSerialNumber() `virtual String ^ GetSerialNumber () [virtual]`

Query the Serial Number of the device.

Returns

The Serial Number.

11.58.3.33 GetSoftwareKey() `array<BYTE> ^ GetSoftwareKey (`
`unsigned int index)`

11.58.3.34 GetSoftwareKeyString() `[1/2] String ^ GetSoftwareKeyString (`
`SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
uint8_t majorversion)`

11.58.3.35 GetSoftwareKeyString() `[2/2] String ^ GetSoftwareKeyString (`
`uint8_t ProgrammID,
uint8_t majorversion)`

11.58.3.36 GetStatus() `virtual uint32_t GetStatus (`
`[System::Runtime::InteropServices::Out] uint32_t% iStatus) [virtual]`

11.58.3.37 GetStatusOfLastCommand() `uint32_t GetStatusOfLastCommand ()`

Gets the status of the last call to the McsUsb Library.

Returns

The Error Status of the last McsUsb command. 0 on success.

11.58.3.38 GetUsbListEntry() `CMcsUsbListEntryNet ^ GetUsbListEntry ()`**11.58.3.39 GetVersion()** [1/2] `DriverVersionNet ^ GetVersion ()`**11.58.3.40 GetVersion()** [2/2] `DriverVersionNet ^ GetVersion (CFirmwareDestinationNet dest)`**11.58.3.41 HasSoftwareKey()** [1/2] `bool HasSoftwareKey (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)`**11.58.3.42 HasSoftwareKey()** [2/2] `bool HasSoftwareKey (uint8_t ProgrammID, uint8_t majorversion)`**11.58.3.43 IsConnected()** `virtual bool IsConnected () [virtual]`

Check if a device is Connected.

Returns

true if the device is connected.

11.58.3.44 IsDeviceHighSpeed() `bool IsDeviceHighSpeed ()`

11.58.3.45 IsDeviceHighSpeedCapable() `bool IsDeviceHighSpeedCapable ()`

11.58.3.46 IsExceptionsEnabled() `bool IsExceptionsEnabled ()`

11.58.3.47 MultibootGetCypressImageId() `uint32_t MultibootGetCypressImageId (unsigned int sector)`

Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.

Returns

The magic ident code of the image.

11.58.3.48 MultibootGetImageId() `String ^ MultibootGetImageId (unsigned int sector)`

Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.

Returns

The magic ident code of the image.

11.58.3.49 MultibootGetSelectedImage() `uint32_t MultibootGetSelectedImage ()`

Gets sector index of selected FPGA boot image on IFB

Returns

Sector index of image; range: 0..2

11.58.3.50 MultibootSelectImage() `void MultibootSelectImage (unsigned int sector)`

Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.

Returns

Throws exception on error.

11.58.3.51 ReadEepromRegisterPreconfig() [1/3] `bool ReadEepromRegisterPreconfig (`
`uint32_t EEPROMBase,`
`uint32_t DMA_reg,`
`[System::Runtime::InteropServices::Out] uint32_t% DMA_value)`

11.58.3.52 ReadEepromRegisterPreconfig() [2/3] `bool ReadEepromRegisterPreconfig (`
`uint32_t EEPROMBase,`
`uint32_t DMA_reg,`
`[System::Runtime::InteropServices::Out] uint32_t% DMA_value,`
`uint32_t EEPROMSize)`

11.58.3.53 ReadEepromRegisterPreconfig() [3/3] `bool ReadEepromRegisterPreconfig (`
`uint32_t EEPROMBase,`
`uint32_t DMA_reg,`
`[System::Runtime::InteropServices::Out] uint32_t% DMA_value,`
`uint32_t EEPROMSize,`
`uint32_t EepromStartAddress)`

11.58.3.54 ReadRegister() [1/2] `unsigned int ReadRegister (`
`unsigned int reg)`

11.58.3.55 ReadRegister() [2/2] `array<uint32_t> ^ ReadRegister (`
`unsigned int reg,`
`int length)`

11.58.3.56 ReadRegister32() `unsigned int ReadRegister32 (`
`unsigned int adr)`

11.58.3.57 ReadRegisterTimeSlot() `unsigned int ReadRegisterTimeSlot (`
`unsigned int reg,`
`int TimeSlot)`

11.58.3.58 RemoveSoftwareKey() `void RemoveSoftwareKey (`
`unsigned int index)`

11.58.3.59 RescanHeadstage() `void RescanHeadstage (`
`uint32_t headstage)`

rescans and activates a headstage

Parameters

in	<i>headstage</i>	the headstage number (0 or 1)
----	------------------	-------------------------------

11.58.3.60 SetConfiguration() void SetConfiguration (
 uint8_t *config*)

11.58.3.61 SetSoftwareKey() void SetSoftwareKey (
 unsigned int *index*,
 array< BYTE >^ *buffer*)

11.58.3.62 ThrowCUsbExceptionNetOnError() void ThrowCUsbExceptionNetOnError (
 uint32_t *status*)

11.58.3.63 TxnGetSerialNumber() unsigned int TxnGetSerialNumber ()

11.58.3.64 TxnSetSerialNumber() void TxnSetSerialNumber (
 unsigned int *number*)

11.58.3.65 TxnTestMemoryReadAndCheck() unsigned int TxnTestMemoryReadAndCheck (
 unsigned short *index*)

11.58.3.66 TxnTestMemoryWrite() unsigned int TxnTestMemoryWrite (
 unsigned short *index*)

11.58.3.67 ValidKey() [1/2] bool ValidKey (
 String^ *key*,
 [System::Runtime::InteropServices::Out] String^% *serial_number*)

11.58.3.68 ValidKey() [2/2] bool ValidKey (
 String^ key,
 uint8_t ProgrammID,
 uint8_t majorversion,
 [System::Runtime::InteropServices::Out] String^% serial_number)

11.58.3.69 WriteEepromRegisterPreconfig() [1/3] void WriteEepromRegisterPreconfig (
 uint32_t EEPROMBase,
 uint32_t DMA_reg,
 uint32_t DMA_value)

11.58.3.70 WriteEepromRegisterPreconfig() [2/3] void WriteEepromRegisterPreconfig (
 uint32_t EEPROMBase,
 uint32_t DMA_reg,
 uint32_t DMA_value,
 uint32_t EEPROMSize)

11.58.3.71 WriteEepromRegisterPreconfig() [3/3] void WriteEepromRegisterPreconfig (
 uint32_t EEPROMBase,
 uint32_t DMA_reg,
 uint32_t DMA_value,
 uint32_t EEPROMSize,
 uint32_t EepromStartAddress)

11.58.3.72 WriteRegister() [1/2] void WriteRegister (
 unsigned int reg,
 array< unsigned int >^ values)

11.58.3.73 WriteRegister() [2/2] void WriteRegister (
 unsigned int reg,
 unsigned int value)

11.58.3.74 WriteRegister32() void WriteRegister32 (
 unsigned int adr,
 unsigned int value)

11.58.3.75 WriteRegisterArray() `void WriteRegisterArray (`
 `unsigned int reg,`
 `array< unsigned int >^ values)`

11.58.3.76 WriteRegisterTimeSlot() [1/2] `void WriteRegisterTimeSlot (`
 `unsigned int reg,`
 `array< unsigned int >^ values,`
 `int TimeSlot)`

11.58.3.77 WriteRegisterTimeSlot() [2/2] `void WriteRegisterTimeSlot (`
 `unsigned int reg,`
 `unsigned int value,`
 `int TimeSlot)`

11.58.3.78 WriteRegisterValue() `void WriteRegisterValue (`
 `unsigned int reg,`
 `unsigned int value)`

11.58.4 Member Data Documentation

11.58.4.1 Status_AlreadyConfigured `const uint32_t Status_AlreadyConfigured = (0xE0110001L)`
[static]

11.58.4.2 Status_BadStartFrame `const uint32_t Status_BadStartFrame = (0xE0100A00L)` [static]

11.58.4.3 Status_Btstuff `const uint32_t Status_Btstuff = (0xE0100002L)` [static]

11.58.4.4 Status_BufferOverrun `const uint32_t Status_BufferOverrun = (0xE010000CL)` [static]

11.58.4.5 Status_BufferUnderrun `const uint32_t Status_BufferUnderrun = (0xE010000DL)` [static]

11.58.4.6 Status_Canceled `const uint32_t Status_Canceled = (0xE0110000L) [static]`

11.58.4.7 Status_Canceling `const uint32_t Status_Canceling = (0xE0120000L) [static]`

11.58.4.8 Status_ConnectedPipes `const uint32_t Status_ConnectedPipes = (0xE01F000AL) [static]`

11.58.4.9 Status_ControlNotOwned `const uint32_t Status_ControlNotOwned = (0xE0100D00L) [static]`

11.58.4.10 Status_Crc `const uint32_t Status_Crc = (0xE0100001L) [static]`

11.58.4.11 Status_DataOverrun `const uint32_t Status_DataOverrun = (0xE0100008L) [static]`

11.58.4.12 Status_DataToggleMismatch `const uint32_t Status_DataToggleMismatch = (0xE0100003L) [static]`

11.58.4.13 Status_DataUnderrun `const uint32_t Status_DataUnderrun = (0xE0100009L) [static]`

11.58.4.14 Status_DeviceLocked `const uint32_t Status_DeviceLocked = (0xE01F0010L) [static]`

11.58.4.15 Status_DeviceNotFound `const uint32_t Status_DeviceNotFound = (0xE01F0003L) [static]`

11.58.4.16 Status_DeviceRemoved `const uint32_t Status_DeviceRemoved = (0xE01F0008L) [static]`

11.58.4.17 Status_DevNotResponding `const uint32_t Status_DevNotResponding = (0xE0100005L)`
[static]

11.58.4.18 Status_EndpointHalted `const uint32_t Status_EndpointHalted = (0xE0100030L)` [static]

11.58.4.19 Status_ErrorBusy `const uint32_t Status_ErrorBusy = (0xE0100400L)` [static]

11.58.4.20 Status_ErrorShortTransfer `const uint32_t Status_ErrorShortTransfer = (0xE0100900L)`
[static]

11.58.4.21 Status_Fifo `const uint32_t Status_Fifo = (0xE0100010L)` [static]

11.58.4.22 Status_FrameControlOwned `const uint32_t Status_FrameControlOwned = (0xE0100C00L)`
[static]

11.58.4.23 Status_InternalHcError `const uint32_t Status_InternalHcError = (0xE0100800L)` [static]

11.58.4.24 Status_InvalidDeviceHandle `const uint32_t Status_InvalidDeviceHandle = (0xE0100013L)`
[static]

11.58.4.25 Status_InvalidHandle `const uint32_t Status_InvalidHandle = (0xE0100012L)` [static]

11.58.4.26 Status_InvalidParameter `const uint32_t Status_InvalidParameter = (0xE0100300L)` [static]

11.58.4.27 Status_InvalidPipeHandle `const uint32_t Status_InvalidPipeHandle = (0xE0100600L)`
[static]

11.58.4.28 Status_InvalidUrbFunction `const uint32_t Status_InvalidUrbFunction = (0xE0100200L)`
[static]

11.58.4.29 Status_IoPending `const uint32_t Status_IoPending = (0xE01F0006L)` [static]

11.58.4.30 Status_IoTimeout `const uint32_t Status_IoTimeout = (0xE01F0007L)` [static]

11.58.4.31 Status_IsochRequestFailed `const uint32_t Status_IsochRequestFailed = (0xE0100B00L)`
[static]

11.58.4.32 Status_LastUsbErrorMismatch `const uint32_t Status_LastUsbErrorMismatch = (0xE01F0022L)` [static]

11.58.4.33 Status_NoBandwidth `const uint32_t Status_NoBandwidth = (0xE0100700L)` [static]

11.58.4.34 Status_NoMemory `const uint32_t Status_NoMemory = (0xE0100100L)` [static]

11.58.4.35 Status_NoSuchDevice `const uint32_t Status_NoSuchDevice = (0xE01F0002L)` [static]

11.58.4.36 Status_NotAccessed `const uint32_t Status_NotAccessed = (0xE010000FL)` [static]

11.58.4.37 Status_NotSupported `const uint32_t Status_NotSupported = (0xE01F0005L)` [static]

11.58.4.38 Status_PidCheckFailure `const uint32_t Status_PidCheckFailure = (0xE0100006L) [static]`

11.58.4.39 Status_PipeNotLinked `const uint32_t Status_PipeNotLinked = (0xE01F0009L) [static]`

11.58.4.40 Status_RequestFailed `const uint32_t Status_RequestFailed = (0xE0100500L) [static]`

11.58.4.41 Status_RequestMutexFailed `const uint32_t Status_RequestMutexFailed = (0xE01F0021L) [static]`

11.58.4.42 Status_RequestMutexTimeout `const uint32_t Status_RequestMutexTimeout = (0xE01↵F0020L) [static]`

11.58.4.43 Status_Stall `const uint32_t Status_Stall = (0xE0100004L) [static]`

11.58.4.44 Status_Unconfigured `const uint32_t Status_Unconfigured = (0xE0110002L) [static]`

11.58.4.45 Status_UnexpectedPid `const uint32_t Status_UnexpectedPid = (0xE0100007L) [static]`

11.58.4.46 WPAError_ScanningIsPending `const uint32_t WPAError_ScanningIsPending = ((0x↵A0220000L) | 0x0036) [static]`


11.58.5 Property Documentation

11.58.5.1 SerialNumber `virtual String^ SerialNumber [get]`

11.59 CMcsUsbPointerContainer Class Reference

11.60 CMEA2100_256DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CMEA2100_256DacqGroupChannelSelectionNet:

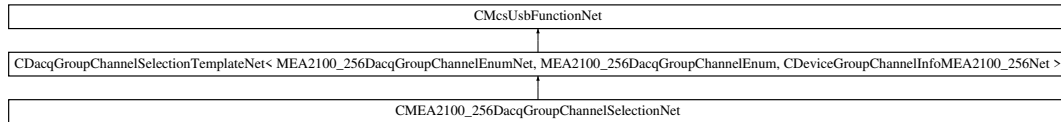

Public Member Functions

- [CMEA2100_256DacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)

Additional Inherited Members


11.60.1 Constructor & Destructor Documentation

11.60.1.1 CMEA2100_256DacqGroupChannelSelectionNet() [CMEA2100_256DacqGroupChannelSelectionNet](#)
(
 [CMcsUsbNet](#)[^] mcsusb)

11.61 CMEA2100x256FunctionNet Class Reference

[CMEA2100x256FunctionNet](#) is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

Inheritance diagram for CMEA2100x256FunctionNet:

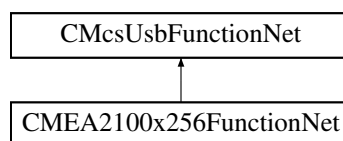

Public Member Functions

- [CMEA2100x256FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pMEA2100x256↔
FunctionPointerContainer)

Initializes a new instance of the [CMEA2100x256FunctionNet](#) class.

- [CMEA2100x256FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMEA2100x256FunctionNet](#) ()
- [!CMEA2100x256FunctionNet](#) ()
- [StimulationLayoutConfigurationEnumNet](#) [GetLayoutConfiguration](#) ()

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulat](#)

- void [SetLayoutConfiguration](#) ([StimulationLayoutConfigurationEnumNet](#) LayoutConfiguration)

Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulat](#)

Additional Inherited Members

11.61.1 Detailed Description

[CMEA2100x256FunctionNet](#) is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

11.61.2 Constructor & Destructor Documentation

11.61.2.1 CMEA2100x256FunctionNet() [1/2] [CMEA2100x256FunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*,
 [CMcsUsbFunctionPointerContainer](#)[^] *pMEA2100x256FunctionPointerContainer*)

Initializes a new instance of the [CMEA2100x256FunctionNet](#) class.

11.61.2.2 CMEA2100x256FunctionNet() [2/2] [CMEA2100x256FunctionNet](#) (
 [CMcsUsbNet](#)[^] *mcsusb*)

11.61.2.3 ~CMEA2100x256FunctionNet() [virtual ~CMEA2100x256FunctionNet](#) () [virtual]

11.61.2.4 !"CMEA2100x256FunctionNet() [!CMEA2100x256FunctionNet](#) ()

11.61.3 Member Function Documentation

11.61.3.1 GetLayoutConfiguration() [StimulationLayoutConfigurationEnumNet](#) [GetLayoutConfiguration](#) ()

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).

Returns

The currently active stimulation layout configuration.

11.61.3.2 SetLayoutConfiguration() [void SetLayoutConfiguration](#) (
 [StimulationLayoutConfigurationEnumNet](#) *LayoutConfiguration*)


Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is [Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels](#) divided by [Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode](#).

Parameters

<i>LayoutConfiguration</i>	The new stimulation layout configuration.
----------------------------	---

11.62 CMeaAudioFunctionNet Class Reference

Inheritance diagram for CMeaAudioFunctionNet:

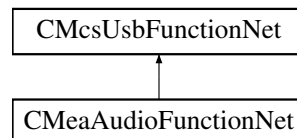

Classes

- struct [s_setaudionet](#)

Public Member Functions

- [CMeaAudioFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meaAudioFunctionPointerContainer)
- [CMeaAudioFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual uint32_t [GetNumberOfAudioChannels](#) ()
Gets the number of available audio channels.
- virtual uint32_t [SetAudioChannels](#) (array< [s_setaudionet](#)[^]>[^] channels)
Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [SetAudioChannels](#) (array< [s_setaudionet](#)[^]>[^] channels, unsigned int virtualDevice)
Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [GetAudioChannels](#) ([System::Runtime::InteropServices::Out]array< [s_setaudionet](#)[^]>[^]% channels)
Gets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t [GetAudioChannels](#) ([System::Runtime::InteropServices::Out]array< [s_setaudionet](#)[^]>[^]% channels, unsigned int virtualDevice)
Gets the electrode to monitor and amplification for the audio channels.

Additional Inherited Members

11.62.1 Constructor & Destructor Documentation

- 11.62.1.1 CMeaAudioFunctionNet()** [1/2] [CMeaAudioFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] meaAudioFunctionPointerContainer)

11.62.1.2 CMeaAudioFunctionNet() [2/2] `CMeaAudioFunctionNet (CMcsUsbNet^ mcsusb)`

11.62.2 Member Function Documentation

11.62.2.1 GetAudioChannels() [1/2] `virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out] array< s_setaudionet^>^% channels) [virtual]`

Gets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which contains the electrode (channel) and amplification on return.
-----------------	--

Returns

Error Status. 0 on success.

11.62.2.2 GetAudioChannels() [2/2] `virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out] array< s_setaudionet^>^% channels, unsigned int virtualDevice) [virtual]`

Gets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which contains the electrode (channel) and amplification on return.
-----------------	--

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------

Returns

Error Status. 0 on success.

11.62.2.3 GetNumberOfAudioChannels() `virtual uint32_t GetNumberOfAudioChannels () [virtual]`

Gets the number of available audio channels.

Returns

The number of audio channels available, 0 when there are none.

11.62.2.4 SetAudioChannels() `[1/2] virtual uint32_t SetAudioChannels (array< s_setaudionet^>^ channels) [virtual]`

Sets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which defines the electrode (channel) and amplification.
-----------------	---

Returns

Error Status. 0 on success.

11.62.2.5 SetAudioChannels() `[2/2] virtual uint32_t SetAudioChannels (array< s_setaudionet^>^ channels, unsigned int virtualDevice) [virtual]`

Sets the electrode to monitor and amplification for the audio channels.

Parameters

<i>channels</i>	Struct which defines the electrode (channel) and amplification.
-----------------	---

Parameters

<i>virtualDevice</i>	Virtual device to use.
----------------------	------------------------


Returns

Error Status. 0 on success.

11.63 CMeaCleanDeviceNet Class Reference

[CMeaCleanDeviceNet](#) is the class to access the MEA Clean device.

Inheritance diagram for CMeaCleanDeviceNet:

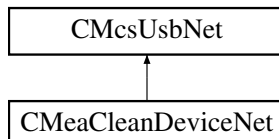

Public Member Functions

- [CMeaCleanDeviceNet](#) ()
Initializes a new instance of the [CMeaCleanDeviceNet](#) class.
- virtual [~CMeaCleanDeviceNet](#) ()
- [!CMeaCleanDeviceNet](#) ()
- void [Start](#) ()
Starts a MEA Clean run.
- void [Stop](#) ()
Stops a MEA Clean run.
- void [SetSlope](#) (uint32_t voltageSlope)
Sets the voltage slope.
- void [SetCycles](#) (uint32_t cycles)
Sets the number of cycles.
- void [SetMinVoltage](#) (int32_t voltageMin)
Sets the lower voltage level.
- void [SetMaxVoltage](#) (int32_t voltageMax)
Sets the upper voltage level.
- bool [IsRunning](#) ()
Gets if the MEA Clean device is running.
- uint32_t [GetSlope](#) ()
Gets the voltage slope.
- uint32_t [GetCycles](#) ()
Gets the number of cycles.
- int32_t [GetMinVoltage](#) ()
Gets the lower voltage level.
- int32_t [GetMaxVoltage](#) ()
Gets the upper voltage level.
- int32_t [GetOutputVoltage](#) ()
Gets the output voltage.
- int32_t [GetCycle](#) ()
Gets the current cycle.

Additional Inherited Members

11.63.1 Detailed Description

[CMeaCleanDeviceNet](#) is the class to access the MEA Clean device.

11.63.2 Constructor & Destructor Documentation

11.63.2.1 CMeaCleanDeviceNet() [CMeaCleanDeviceNet](#) ()

Initializes a new instance of the [CMeaCleanDeviceNet](#) class.

11.63.2.2 ~CMeaCleanDeviceNet() [virtual ~CMeaCleanDeviceNet](#) () [virtual]

11.63.2.3 "!CMeaCleanDeviceNet() [!CMeaCleanDeviceNet](#) ()

11.63.3 Member Function Documentation

11.63.3.1 GetCycle() [int32_t GetCycle](#) ()

Gets the current cycle.

Returns

The cycle number.

11.63.3.2 GetCycles() [uint32_t GetCycles](#) ()

Gets the number of cycles.

Returns

The number of cycles to run for.

11.63.3.3 GetMaxVoltage() `int32_t GetMaxVoltage ()`

Gets the upper voltage level

Returns

The upper voltage level in mV.

11.63.3.4 GetMinVoltage() `int32_t GetMinVoltage ()`

Gets the lower voltage level.

Returns

The lower voltage level in mV.

11.63.3.5 GetOutputVoltage() `int32_t GetOutputVoltage ()`

Gets the output voltage.

Returns

The output voltage in mV.

11.63.3.6 GetSlope() `uint32_t GetSlope ()`

Gets the voltage slope.

Returns

The voltage slope in mV/s.

11.63.3.7 IsRunning() `bool IsRunning ()`

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

11.63.3.8 SetCycles() `void SetCycles (
 uint32_t cycles)`

Sets the number of cycles.

Parameters

<i>cycles</i>	The number of cycles to run for (0 .. 99).
---------------	--

11.63.3.9 SetMaxVoltage() `void SetMaxVoltage (`
`int32_t voltageMax)`

Sets the upper voltage level.

Parameters

<i>voltageMax</i>	The upper voltage level in mV (-1.6 .. 1.6 V).
-------------------	--

11.63.3.10 SetMinVoltage() `void SetMinVoltage (`
`int32_t voltageMin)`

Sets the lower voltage level.

Parameters

<i>voltageMin</i>	The lower voltage level in mV (-1.6 .. 1.6 V).
-------------------	--

11.63.3.11 SetSlope() `void SetSlope (`
`uint32_t voltageSlope)`

Sets the voltage slope.

Parameters

<i>voltageSlope</i>	The voltage slope in mV/s (range 0 .. 60 V/s).
---------------------	--

11.63.3.12 Start() `void Start ()`

Starts a MEA Clean run.


11.63.3.13 Stop() `void Stop ()`

Stops a MEA Clean run.

11.64 CMeaCoatDeviceNet Class Reference

[CMeaCoatDeviceNet](#) is the class to access the MEA Coat device.

Inheritance diagram for CMeaCoatDeviceNet:

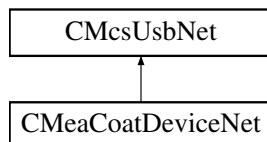

Public Member Functions

- [CMeaCoatDeviceNet](#) ()
Initializes a new instance of the [CMeaCoatDeviceNet](#) class.
- virtual [~CMeaCoatDeviceNet](#) ()
- [!CMeaCoatDeviceNet](#) ()
- void [Start](#) ()
Starts a MEA Coat run.
- void [Stop](#) ()
Stops a MEA Coat run.
- void [SetSlope](#) (int32_t currentSlope)
Sets the current slope.
- void [SetDuration](#) (uint32_t duration)
Sets the duration of a MEA Coat run.
- void [SetMaxCurrent](#) (uint32_t currentMax)
Sets the limit of the current ramp (absolute value).
- void [SetOffsetCurrent](#) (int32_t currentOffset)
Sets the offset of the current.
- bool [IsRunning](#) ()
Gets if the MEA Clean device is running.
- int32_t [GetSlope](#) ()
Gets the current slope.
- uint32_t [GetDuration](#) ()
Gets the duration of a MEA Coat run.
- uint32_t [GetMaxCurrent](#) ()
Gets the limit of the current ramp (absolute value).
- int32_t [GetOffsetCurrent](#) ()
Gets the offset of the current.
- int32_t [GetOutputCurrent](#) ()
Gets the output current.
- int32_t [GetTimeInPlateau](#) ()
Gets the time in the plateau.
- void [SetPauseDuration](#) (uint32_t pauseDuration)
Sets the duration of the pause between MEA Coat pulses.
- uint32_t [GetPauseDuration](#) ()
Gets the duration of the pause between MEA Coat pulses.
- int32_t [GetTimeInPause](#) ()
Gets the time in the pause.

- void [SetCycles](#) (uint32_t cycles)
Sets the number of cycles.
- uint32_t [GetCycles](#) ()
Gets the number of cycles.
- int32_t [GetCurrentCycle](#) ()
Gets the current cycle.

Additional Inherited Members

11.64.1 Detailed Description

[CMeaCoatDeviceNet](#) is the class to access the MEA Coat device.

11.64.2 Constructor & Destructor Documentation

11.64.2.1 CMeaCoatDeviceNet() [CMeaCoatDeviceNet](#) ()

Initializes a new instance of the [CMeaCoatDeviceNet](#) class.

11.64.2.2 ~CMeaCoatDeviceNet() [virtual](#) [~CMeaCoatDeviceNet](#) () [\[virtual\]](#)

11.64.2.3 "!CMeaCoatDeviceNet() [!CMeaCoatDeviceNet](#) ()

11.64.3 Member Function Documentation

11.64.3.1 GetCurrentCycle() [int32_t](#) [GetCurrentCycle](#) ()

Gets the current cycle.

Returns

The cycle number.

11.64.3.2 GetCycles() `uint32_t GetCycles ()`

Gets the number of cycles.

Returns

The number of cycles to run for.

11.64.3.3 GetDuration() `uint32_t GetDuration ()`

Gets the duration of a MEA Coat run.

Returns

The duration in ms.

11.64.3.4 GetMaxCurrent() `uint32_t GetMaxCurrent ()`

Gets the limit of the current ramp (absolute value).

Returns

The limit of the current ramp in pA (absolute value).

11.64.3.5 GetOffsetCurrent() `int32_t GetOffsetCurrent ()`

Gets the offset of the current.

Returns

The offset of the current in pA.

11.64.3.6 GetOutputCurrent() `int32_t GetOutputCurrent ()`

Gets the output current.

Returns

The output current in pA.

11.64.3.7 GetPauseDuration() `uint32_t GetPauseDuration ()`

Gets the duration of the pause between MEA Coat pulses.

Returns

The duration in ms.

11.64.3.8 GetSlope() `int32_t GetSlope ()`

Gets the current slope.

Returns

The current slope in pA/s.

11.64.3.9 GetTimeInPause() `int32_t GetTimeInPause ()`

Gets the time in the pause.

Returns

The time in the pause in ms.

11.64.3.10 GetTimeInPlateau() `int32_t GetTimeInPlateau ()`

Gets the time in the plateau.

Returns

The time in the plateau in ms.

11.64.3.11 IsRunning() `bool IsRunning ()`

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

11.64.3.12 SetCycles() `void SetCycles (
uint32_t cycles)`

Sets the number of cycles.

Parameters

<i>cycles</i>	The number of cycles to run for (0 .. 99).
---------------	--

11.64.3.13 SetDuration() `void SetDuration (`
 `uint32_t duration)`

Sets the duration of a MEA Coat run.

Parameters

<i>duration</i>	The duration in ms (range 0 .. 65 s).
-----------------	---------------------------------------

11.64.3.14 SetMaxCurrent() `void SetMaxCurrent (`
 `uint32_t currentMax)`

Sets the limit of the current ramp (absolute value).

Parameters

<i>currentMax</i>	The limit of the current ramp in pA (absolute value, 0 .. 18 nA).
-------------------	---

11.64.3.15 SetOffsetCurrent() `void SetOffsetCurrent (`
 `int32_t currentOffset)`

Sets the offset of the current.

Parameters

<i>currentOffset</i>	The offset of the current in pA (-10 .. 10 nA).
----------------------	---

11.64.3.16 SetPauseDuration() `void SetPauseDuration (`
 `uint32_t pauseDuration)`

Sets the duration of the pause between MEA Coat pulses.

Parameters

<i>pauseDuration</i>	The duration in ms (range 0 .. 65 s).
----------------------	---------------------------------------

11.64.3.17 SetSlope() `void SetSlope (`
`int32_t currentSlope)`

Sets the current slope.

Parameters

<i>currentSlope</i>	The current slope in pA/s (range -65 .. 65 nA/s).
---------------------	---

11.64.3.18 Start() `void Start ()`

Starts a MEA Coat run.


11.64.3.19 Stop() `void Stop ()`

Stops a MEA Coat run.

11.65 CMeaDeviceNet Class Reference

Base class for MEA data acquisition devices.

Inheritance diagram for CMeaDeviceNet:

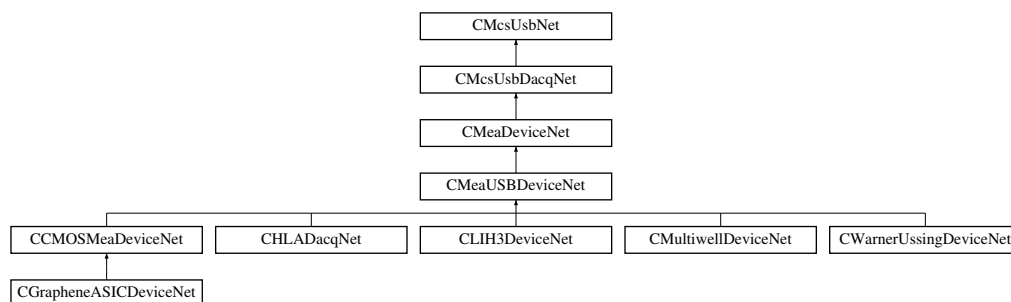

Public Member Functions

- [CMeaDeviceNet](#) ([McsBusTypeEnumNet](#) bustype)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [CMeaDeviceNet](#) ([McsBusTypeEnumNet](#) bustype, [OnChannelData](#)[^] channelData, [OnError](#)[^] error)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [~CMeaDeviceNet](#) ()
- virtual `int32_t` [GetGain](#) ()
Gets the amplifier gain of the device.

- `int32_t GetEnumerationSpeed ()`
- `virtual int32_t GetAnalogGain ()`
Gets the gain of the analog inputs of the device.
- `virtual uint32_t EnableDigitalIn (bool enable, unsigned int virtualDevice)`
Enable the digital data word in the datastream.
- `virtual uint32_t EnableDigitalIn (DigitalDatastreamEnableEnumNet enable, unsigned int virtualDevice)`
Enable digital data words in the datastream.
- `virtual uint32_t EnableTimestamp (bool enable, unsigned int virtualDevice)`
Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.
- `virtual uint32_t EnableChecksum (bool enable, unsigned int virtualDevice)`
Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.
- `virtual void SetDigitalOut (unsigned int digout_value, int pulselength)`
Generate a pulse on the digital output.
- `virtual uint32_t SetNumberOfChannels (int NumberOfChannels)`
Sets the number of analog channels in the datastream.
- `virtual uint32_t SetNumberOfChannels (int NumberOfChannels, unsigned int virtualDevice)`
Sets the number of analog channels in the datastream.
- `virtual uint32_t SetNumberOfAnalogChannels (unsigned int NumberOfChannels_HS1, unsigned int NumberOfChannels_HS2, unsigned int NumberOfChannels_DSP, unsigned int NumberOfChannels_IF, unsigned int virtualDevice)`
Sets the number of analog channels in the datastream for the MEA2100 device.
- `virtual uint32_t SetTriggerPeriod (int samples, unsigned int virtualDevice)`
Sets the maximum number of samples per trigger.
- `virtual uint32_t SetTriggerMaskValue (unsigned int mask, unsigned int value, unsigned int virtualDevice)`
Defines a pattern on the digital dataword which will start a trigger when found.

Properties

- `CMeFunctionNet^ MeFunctionNet [get]`
- `CWClassicFunctionNet^ WClassicFunctionNet [get]`
- `CW2100_FunctionNet^ W2100_FunctionNet [get]`
- `CMeaAudioFunctionNet^ MeaAudioFunctionNet [get]`
- `CMeaDigitalDataFunctionNet^ MeaDigitalDataFunctionNet [get]`
- `CMeaFeedbackFunctionNet^ MeaFeedbackFunctionNet [get]`
- `virtual int Gain [get]`
The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.
- `virtual int AnalogGain [get]`
The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Additional Inherited Members

11.65.1 Detailed Description

Base class for MEA data acquisition devices.

There are two different device types for MEA data acquisition devices. There are the USB-MEA devices and the MC↔_Card. In .NET both classes can be accessed by the constructor of the base class `CMeaDeviceNet`, which constructs the correct underlying C++ class for the USB-MEA device on the one hand or the MC_Card device on the other hand. Through this interface both device types USB-MEA devices and MC_Card devices can be accessed

11.65.2 Constructor & Destructor Documentation

11.65.2.1 CMeaDeviceNet() [1/2] CMeaDeviceNet (McsBusTypeEnumNet bustype)

Initializes a new instance of CMeaDeviceNet class.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

11.65.2.2 CMeaDeviceNet() [2/2] CMeaDeviceNet (McsBusTypeEnumNet bustype, OnChannelData^ channelData, OnError^ error)

Initializes a new instance of CMeaDeviceNet class.

Parameters

<i>bustype</i>	Type of device to use, either USB or PCI.
----------------	---

Parameters

<i>channelData</i>	Callback to call when new data is available.
--------------------	--

Parameters

<i>error</i>	Callback to call when an error occurred.
--------------	--

11.65.2.3 ~CMeaDeviceNet() ~CMeaDeviceNet ()

11.65.3 Member Function Documentation

11.65.3.1 EnableChecksum() `virtual uint32_t EnableChecksum (`
 `bool enable,`
 `unsigned int virtualDevice) [virtual]`

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

11.65.3.2 EnableDigitalIn() `[1/2] virtual uint32_t EnableDigitalIn (`
 `bool enable,`
 `unsigned int virtualDevice) [virtual]`

Enable the digital data word in the datastream.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

11.65.3.3 EnableDigitalIn() `[2/2] virtual uint32_t EnableDigitalIn (`
 `DigitalDatastreamEnableEnumNet enable,`
 `unsigned int virtualDevice) [virtual]`

Enable digital data words in the datastream.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

11.65.3.4 EnableTimestamp() `virtual uint32_t EnableTimestamp (`
 `bool enable,`
 `unsigned int virtualDevice) [virtual]`

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

Parameters

<i>enable</i>	True to enable, False to disable.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

11.65.3.5 GetAnalogGain() `virtual int32_t GetAnalogGain () [virtual]`

Gets the gain of the analog inputs of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.65.3.6 GetEnumerationSpeed() `int32_t GetEnumerationSpeed ()`

11.65.3.7 GetGain() `virtual int32_t GetGain () [virtual]`

Gets the amplifier gain of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.65.3.8 SetDigitalOut() `virtual void SetDigitalOut (`
 `unsigned int digout_value,`
 `int pulselength) [virtual]`

Generate a pulse on the digital output.

Parameters

<i>digout_value</i>	Bitmask to set on the digital out.
---------------------	------------------------------------

Parameters

<i>pulselength</i>	Pulselength in ms.
--------------------	--------------------

11.65.3.9 SetNumberOfAnalogChannels() `virtual uint32_t SetNumberOfAnalogChannels (`
 `unsigned int NumberOfChannels_HS1,`
 `unsigned int NumberOfChannels_HS2,`
 `unsigned int NumberOfChannels_DSP,`
 `unsigned int NumberOfChannels_IF,`
 `unsigned int virtualDevice) [virtual]`

Sets the number of analog channels in the datastream for the MEA2100 device.

Parameters

<i>NumberOfChannels_HS1</i>	Number of analog channels from the Headstage 1.
-----------------------------	---

Parameters

<i>NumberOfChannels_HS2</i>	Number of analog channels from the Headstage 2.
-----------------------------	---

Parameters

<i>NumberOfChannels_DSP</i>	Number of data words from the DSP.
-----------------------------	------------------------------------

Parameters

<i>NumberOfChannels↵_IF</i>	Number of analog channels from the Interfaceboard.
-----------------------------	--

Parameters

<i>virtualDevice</i>	virtualDevice to use.
----------------------	-----------------------

Returns

Error Status. 0 on success.

11.65.3.10 SetNumberOfChannels() [1/2] virtual uint32_t SetNumberOfChannels (
int *NumberOfChannels*) [virtual]

Sets the number of analog channels in the datastream.

Parameters

<i>NumberOfChannels</i>	Number of analog channels.
-------------------------	----------------------------

Returns

Error Status. 0 on success.

11.65.3.11 SetNumberOfChannels() [2/2] virtual uint32_t SetNumberOfChannels (
int *NumberOfChannels*,
unsigned int *virtualDevice*) [virtual]

Sets the number of analog channels in the datastream.

Parameters

<i>NumberOfChannels</i>	Number of analog channels.
<i>virtualDevice</i>	virtual device to use.

Returns

Error Status. 0 on success.

11.65.3.12 SetTriggerMaskValue() `virtual uint32_t SetTriggerMaskValue (`
 `unsigned int mask,`
 `unsigned int value,`
 `unsigned int virtualDevice) [virtual]`

Defines a pattern on the digital dataword which will start a trigger when found.

Parameters

<i>mask</i>	Bits in the digital dataword which are monitored for a match with value.
-------------	--

Parameters

<i>value</i>	Pattern which must match for the trigger to start.
--------------	--

Returns

Error Status. 0 on success.

11.65.3.13 SetTriggerPeriod() `virtual uint32_t SetTriggerPeriod (`
 `int samples,`
 `unsigned int virtualDevice) [virtual]`

Sets the maximum number of samples per trigger.

Parameters

<i>samples</i>	Number of samples to acquire after the trigger condition is met.
----------------	--

Returns

Error Status. 0 on success.

11.65.4 Property Documentation

11.65.4.1 AnalogGain `virtual int AnalogGain [get]`


The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.65.4.2 Gain `virtual int Gain [get]`

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.65.4.3 MeaAudioFunctionNet `CMeaAudioFunctionNet^ MeaAudioFunctionNet [get]`**11.65.4.4 MeaDigitalDataFunctionNet** `CMeaDigitalDataFunctionNet^ MeaDigitalDataFunctionNet [get]`**11.65.4.5 MeaFeedbackFunctionNet** `CMeaFeedbackFunctionNet^ MeaFeedbackFunctionNet [get]`**11.65.4.6 MeFunctionNet** `CMeFunctionNet^ MeFunctionNet [get]`**11.65.4.7 W2100_FunctionNet** `CW2100_FunctionNet^ W2100_FunctionNet [get]`**11.65.4.8 WClassicFunctionNet** `CWClassicFunctionNet^ WClassicFunctionNet [get]`**11.66 CMeaDigitalDataFunctionNet Class Reference**

Inheritance diagram for CMeaDigitalDataFunctionNet:

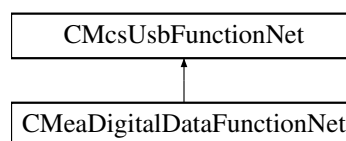**Public Member Functions**

- `CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ meaDigitalFunctionPointerContainer)`
- `CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb)`
- `void SetDigitalData (unsigned int digital_value, unsigned int digital_value_mask)`
Generate a value on the digital output.
- `void SetDigitalData (unsigned int bit_number, bool value)`
Generate a value on the digital output.
- `unsigned int GetDigitalData ()`
Get the value of the digital output.

Additional Inherited Members

11.66.1 Constructor & Destructor Documentation

11.66.1.1 CMeaDigitalDataFunctionNet() [1/2] `CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ meaDigitalFunctionPointerContainer)`

11.66.1.2 CMeaDigitalDataFunctionNet() [2/2] `CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb)`

11.66.2 Member Function Documentation

11.66.2.1 GetDigitalData() `unsigned int GetDigitalData ()`

Get the value of the digital output.

Returns

Value on the digital data register.

11.66.2.2 SetDigitalData() [1/2] `void SetDigitalData (unsigned int bit_number, bool value)`

Generate a value on the digital output.

Parameters

<i>bit_number</i>	Bit number to change.
-------------------	-----------------------

Parameters

<i>value</i>	Bit value.
--------------	------------

11.66.2.3 SetDigitalData() [2/2] void SetDigitalData (
 unsigned int *digital_value*,
 unsigned int *digital_value_mask*)

Generate a value on the digital output.

Parameters


<i>digital_value</i>	Value to set.
----------------------	---------------

Parameters

<i>digital_value_mask</i>	Mask for change.
---------------------------	------------------

11.67 CMeaFeedbackFunctionNet Class Reference

Inheritance diagram for CMeaFeedbackFunctionNet:

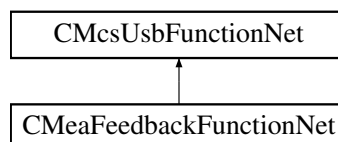

Public Member Functions

- [CMeaFeedbackFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meaFeedback↔
FunctionNet)
- [CMeaFeedbackFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [FeedbackSetFeedback](#) (unsigned char on, unsigned short digoutmask, unsigned short diginmask)
- unsigned int [FeedbackGetSampleTimerCount](#) ([System::Runtime::InteropServices::Out]unsigned int% CurrentCount, [System::Runtime::InteropServices::Out]unsigned int% LastKnownCount, [System::Runtime↔
::InteropServices::Out]bool% On)
- void [FeedbackSetDigitalMapping](#) (unsigned short channel, unsigned short outmapping, unsigned short in-
mapping)
- void [FeedbackSetFilterParameter](#) (unsigned char filter, array< short >[^] parameters)
- void [FeedbackSetFilterParameter32](#) (unsigned char filter, array< int >[^] parameters)
- void [FeedbackSetIIRFilterParameter](#) (unsigned char filter, int length, array< double >[^] parameters)
- void [FeedbackSetMkFilter](#) (unsigned char filter, String[^] filtertype, double cheb_ripple, String[^] passtype, int
order, double alpha1, double alpha2)
- void [FeedbackSetChannelFilter](#) (short channel, char filter)
- void [FeedbackSetGlobalChannelFilter](#) (char filter, unsigned short firstchannel, unsigned short lastchannel)

- void [FeedbackSetFilterOff](#) ()
- void [FeedbackSetNumberOfSpikeDetectors](#) (unsigned short number)
- void [FeedbackSetSpikeDetectorThreshold](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, int threshold1, int threshold2, short slope)
- void [FeedbackSetNumberOfRateCounter](#) (unsigned short number)
- void [FeedbackSetRateCounter](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel)
- void [FeedbackSetNumberOfRateDetectors](#) (unsigned short number)
- void [FeedbackSetRateDetector](#) (unsigned short position, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, unsigned short pulses, unsigned int duration1, unsigned int duration2)
- void [FeedbackSetNumberOfLogics](#) (unsigned short number)
- void [FeedbackSetLogic](#) (unsigned short position, array< unsigned short >^ sourcechannel, unsigned short resultchannel, unsigned int lookup)
- void [FeedbackSetNumberOfTriggers](#) (unsigned short number)
- void [FeedbackSetTrigger](#) (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit)
- void [FeedbackSetAnalogSource](#) ([AnalogSourceEnumNet](#) AnalogSource, unsigned int Channels, unsigned int Offset)

Additional Inherited Members

11.67.1 Constructor & Destructor Documentation

11.67.1.1 CMeaFeedbackFunctionNet() [1/2] [CMeaFeedbackFunctionNet](#) (
[CMcsUsbNet](#)^ *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)^ *meaFeedbackFunctionNet*)

11.67.1.2 CMeaFeedbackFunctionNet() [2/2] [CMeaFeedbackFunctionNet](#) (
[CMcsUsbNet](#)^ *mcsusb*)

11.67.2 Member Function Documentation

11.67.2.1 FeedbackGetSampleTimerCount() unsigned int [FeedbackGetSampleTimerCount](#) (
[System::Runtime::InteropServices::Out] unsigned int% *CurrentCount*,
[System::Runtime::InteropServices::Out] unsigned int% *LastKnownCount*,
[System::Runtime::InteropServices::Out] bool% *On*)

11.67.2.2 FeedbackSetAnalogSource() void [FeedbackSetAnalogSource](#) (
[AnalogSourceEnumNet](#) *AnalogSource*,
unsigned int *Channels*,
unsigned int *Offset*)

11.67.2.3 FeedbackSetChannelFilter() void FeedbackSetChannelFilter (
 short *channel*,
 char *filter*)

11.67.2.4 FeedbackSetDigitalMapping() void FeedbackSetDigitalMapping (
 unsigned short *channel*,
 unsigned short *outmapping*,
 unsigned short *inmapping*)

11.67.2.5 FeedbackSetFeedback() void FeedbackSetFeedback (
 unsigned char *on*,
 unsigned short *digoutmask*,
 unsigned short *diginmask*)

11.67.2.6 FeedbackSetFilterOff() void FeedbackSetFilterOff ()

11.67.2.7 FeedbackSetFilterParameter() void FeedbackSetFilterParameter (
 unsigned char *filter*,
 array< short >^ *parameters*)

11.67.2.8 FeedbackSetFilterParameter32() void FeedbackSetFilterParameter32 (
 unsigned char *filter*,
 array< int >^ *parameters*)

11.67.2.9 FeedbackSetGlobalChannelFilter() void FeedbackSetGlobalChannelFilter (
 char *filter*,
 unsigned short *firstchannel*,
 unsigned short *lastchannel*)

11.67.2.10 FeedbackSetIIRFilterParameter() void FeedbackSetIIRFilterParameter (
 unsigned char *filter*,
 int *length*,
 array< double >^ *parameters*)

11.67.2.11 FeedbackSetLogic() void FeedbackSetLogic (
 unsigned short *position*,
 array< unsigned short >^ *sourcechannel*,
 unsigned short *resultchannel*,
 unsigned int *lookup*)

11.67.2.12 FeedbackSetMkFilter() void FeedbackSetMkFilter (
 unsigned char *filter*,
 String^ *filtertype*,
 double *cheb_ripple*,
 String^ *passtype*,
 int *order*,
 double *alpha1*,
 double *alpha2*)

11.67.2.13 FeedbackSetNumberOfLogics() void FeedbackSetNumberOfLogics (
 unsigned short *number*)

11.67.2.14 FeedbackSetNumberOfRateCounter() void FeedbackSetNumberOfRateCounter (
 unsigned short *number*)

11.67.2.15 FeedbackSetNumberOfRateDetectors() void FeedbackSetNumberOfRateDetectors (
 unsigned short *number*)

11.67.2.16 FeedbackSetNumberOfSpikeDetectors() void FeedbackSetNumberOfSpikeDetectors (
 unsigned short *number*)

11.67.2.17 FeedbackSetNumberOfTriggers() void FeedbackSetNumberOfTriggers (
 unsigned short *number*)

11.67.2.18 FeedbackSetRateCounter() void FeedbackSetRateCounter (
 unsigned short *position*,
 unsigned short *sourcechannel*,
 unsigned short *resultchannel*)

11.67.2.19 FeedbackSetRateDetector() void FeedbackSetRateDetector (

```

    unsigned short position,
    unsigned short resultchannel,
    unsigned short trigger,
    unsigned short totzeit,
    unsigned short pulses,
    unsigned int duration1,
    unsigned int duration2 )
```

11.67.2.20 FeedbackSetSpikeDetectorThreshold() void FeedbackSetSpikeDetectorThreshold (

```

    unsigned short position,
    unsigned short sourcechannel,
    unsigned short resultchannel,
    unsigned short trigger,
    unsigned short totzeit,
    int threshold1,
    int threshold2,
    short slope )
```


11.67.2.21 FeedbackSetTrigger() void FeedbackSetTrigger (

```

    unsigned short position,
    unsigned short sourcechannel,
    unsigned short resultchannel,
    unsigned short trigger,
    unsigned short totzeit )
```

11.68 CMealImpedanceDeviceNet Class Reference

Inheritance diagram for CMealImpedanceDeviceNet:

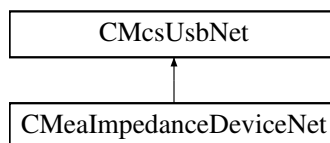

Public Member Functions

- [CMealImpedanceDeviceNet](#) ()
- [~CMealImpedanceDeviceNet](#) ()
- virtual void [StartMeasurement](#) (unsigned short channel)
- virtual unsigned short [GetReady](#) ()
- virtual unsigned short [GetArraySize](#) ()
- virtual array< unsigned short > ^ [GetResult](#) ()
- unsigned short [GetAdapterCode](#) ()
- virtual unsigned int [GetImpedanceTestFrequency](#) ()
- virtual void [SetImpedanceTestFrequency](#) (unsigned int TestFrequency_Hertz)

Additional Inherited Members

11.68.1 Constructor & Destructor Documentation

11.68.1.1 CMeaImpedanceDeviceNet() `CMeaImpedanceDeviceNet ()`

11.68.1.2 ~CMeaImpedanceDeviceNet() `~CMeaImpedanceDeviceNet ()`

11.68.2 Member Function Documentation

11.68.2.1 GetAdapterCode() `unsigned short GetAdapterCode ()`

11.68.2.2 GetArraySize() `virtual unsigned short GetArraySize () [virtual]`

11.68.2.3 GetImpedanceTestFrequency() `virtual unsigned int GetImpedanceTestFrequency () [virtual]`

11.68.2.4 GetReady() `virtual unsigned short GetReady () [virtual]`

11.68.2.5 GetResult() `virtual array<unsigned short> ^ GetResult () [virtual]`


11.68.2.6 SetImpedanceTestFrequency() `virtual void SetImpedanceTestFrequency (unsigned int TestFrequency_Hertz) [virtual]`

11.68.2.7 StartMeasurement() `virtual void StartMeasurement (unsigned short channel) [virtual]`

11.69 CMeasureTableDeviceNet Class Reference

[CMeasureTableDeviceNet](#) is the to control the MCS HLA device

Inheritance diagram for CMeasureTableDeviceNet:

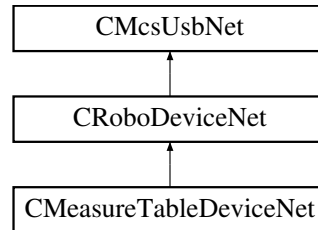

Public Member Functions

- [CMeasureTableDeviceNet](#) (void)

Properties

- [CMcsBus_SensorNet](#)[^] [Sensor](#) [get]

Additional Inherited Members

11.69.1 Detailed Description

[CMeasureTableDeviceNet](#) is the to control the MCS HLA device

11.69.2 Constructor & Destructor Documentation

11.69.2.1 CMeasureTableDeviceNet() [CMeasureTableDeviceNet](#) (
void)


11.69.3 Property Documentation

11.69.3.1 Sensor [CMcsBus_SensorNet](#)[^] [Sensor](#) [get]

11.70 CMeaSwitchDeviceNet Class Reference

The class to control the USB-MEA-Switch.

Inheritance diagram for CMeaSwitchDeviceNet:

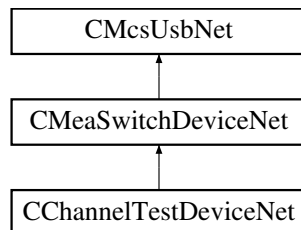

Public Member Functions

- [CMeaSwitchDeviceNet](#) ()
Constructor.
- [~CMeaSwitchDeviceNet](#) ()
Destructor.
- unsigned short [GetNumber](#) ()
Gets the number of boards in the device.
- array< unsigned char > ^ [GetPattern](#) ()
Gets the pattern of the switches that are currently set in the device as char array.
- array< bool > ^ [GetPatternBool](#) ()
Gets the pattern of the switches that are currently set in the device as bools.
- void [SetPattern](#) (array< unsigned char >^ pattern)
Sets the pattern of switches from a char array.
- void [SetPatternBool](#) (array< bool >^ pattern)
Sets the pattern of switches from a.

Additional Inherited Members

11.70.1 Detailed Description

The class to control the USB-MEA-Switch.

This class controls the settings of the USB-MEA-Switch. The box has two inputs for signals from a MEA amplifier. Each of the 64 outputs can be connected to one of the MEAs at the same channel.

11.70.2 Constructor & Destructor Documentation

11.70.2.1 CMeaSwitchDeviceNet() [CMeaSwitchDeviceNet](#) ()

Constructor.

11.70.2.2 `~CMeaSwitchDeviceNet()` `~CMeaSwitchDeviceNet ()`

Destructor.

11.70.3 Member Function Documentation**11.70.3.1** `GetNumber()` `unsigned short GetNumber ()`

Gets the number of boards in the device.

The MEA-Switch are delivered with 64 or 128 channels

11.70.3.2 `GetPattern()` `array<unsigned char> ^ GetPattern ()`

Gets the pattern of the switches that are currently set in the device as char array.

11.70.3.3 `GetPatternBool()` `array<bool> ^ GetPatternBool ()`

Gets the pattern of the switches that are currently set in the device as bools.

11.70.3.4 `SetPattern()` `void SetPattern (array< unsigned char >^ pattern)`

Sets the pattern of switches from a char array.


11.70.3.5 `SetPatternBool()` `void SetPatternBool (array< bool >^ pattern)`

Sets the pattern of switches from a.

11.71 CMeaUSBDeviceNet Class Reference

Class for data acquisition via ME and MEA USB amplifiers

Inheritance diagram for CMeaUSBDeviceNet:

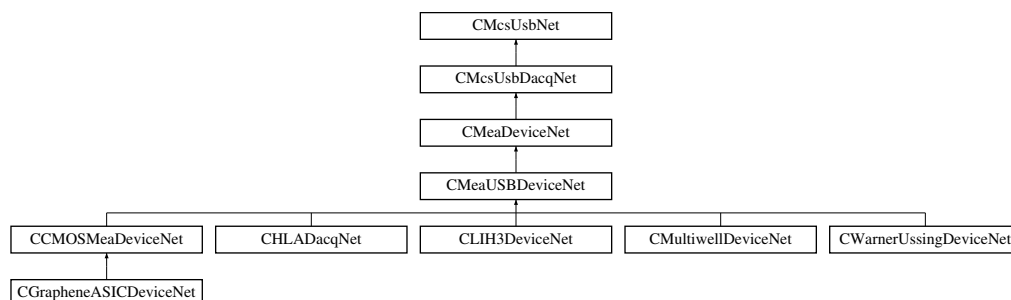

Public Member Functions

- [CMeaUSBDeviceNet](#) ([OnChannelData](#)[^] channelData, [OnError](#)[^] error)
Initializes a new instance of [CMeaDeviceNet](#) class.
- [CMeaUSBDeviceNet](#) ()
Initializes a new instance of [CMeaDeviceNet](#) class.
- [~CMeaUSBDeviceNet](#) ()

Additional Inherited Members

11.71.1 Detailed Description

Class for data acquisition via ME and MEA USB amplifiers

11.71.2 Constructor & Destructor Documentation

11.71.2.1 CMeaUSBDeviceNet() [1/2] [CMeaUSBDeviceNet](#) (
 [OnChannelData](#)[^] channelData,
 [OnError](#)[^] error)

Initializes a new instance of [CMeaDeviceNet](#) class.

Parameters

<i>channelData</i>	Handler to call when new data is available.
--------------------	---

Parameters

<i>error</i>	Handler to call when an error occurs.
--------------	---------------------------------------


11.71.2.2 CMeaUSBDeviceNet() [2/2] [CMeaUSBDeviceNet](#) ()

Initializes a new instance of [CMeaDeviceNet](#) class.

11.71.2.3 ~CMeaUSBDeviceNet() [~CMeaUSBDeviceNet](#) ()

11.72 CMeFunctionNet Class Reference

Inheritance diagram for CMeFunctionNet:

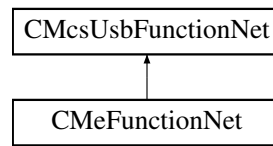

Public Member Functions

- [CMeFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] meFunctionPointer↔ Container)
Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.
- [CMeFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMeFunctionNet](#) (void)
- [!CMeFunctionNet](#) (void)
- void [SetTransformer](#) (unsigned int index, bool onoff)

Additional Inherited Members

11.72.1 Detailed Description

11.72.2 Constructor & Destructor Documentation

11.72.2.1 CMeFunctionNet() [1/2] [CMeFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] meFunctionPointerContainer)

Initializes a new instance of the [CDacCalibrationFunctionNet](#) class.

11.72.2.2 CMeFunctionNet() [2/2] [CMeFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

11.72.2.3 ~CMeFunctionNet() virtual [~CMeFunctionNet](#) (
 void) [virtual]

11.72.2.4 !CMeFunctionNet() [!CMeFunctionNet](#) (
 void)


11.72.3 Member Function Documentation

11.72.3.1 SetTransformer() `void SetTransformer (`
 `unsigned int index,`
 `bool onoff)`

11.73 CMultiBatteryChargerDeviceNet Class Reference

[CMultiBatteryChargerDeviceNet](#) is the class to access the MBC-08 device.

Inheritance diagram for CMultiBatteryChargerDeviceNet:

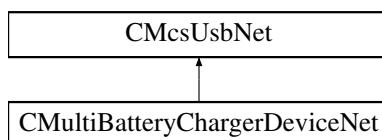

Public Member Functions

- [CMultiBatteryChargerDeviceNet](#) ()
Initializes a new instance of the [CMultiBatteryChargerDeviceNet](#) class.
- virtual [~CMultiBatteryChargerDeviceNet](#) ()
- [!CMultiBatteryChargerDeviceNet](#) ()
- [uint32_t GetChargeCurrent](#) (uint32_t NrChannel)
gets the charge current; unit: mA
- [uint32_t GetDischargeCurrent](#) (uint32_t NrChannel)
gets the discharge current; unit: mA
- void [SetDischargeCurrentSetPoint](#) (uint32_t NrChannel, uint32_t DischargeCurrent_mA)
sets the setpoint for the discharge current; unit: mA
- [uint32_t GetDischargeCurrentSetPoint](#) (uint32_t NrChannel)
gets the setpoint for the discharge current; unit: mA
- void [SetFinalDischargeVoltage](#) (uint32_t NrChannel, uint32_t FinalDischargeVoltage_mV)
sets the final discharge voltage; unit: mV
- [uint32_t GetFinalDischargeVoltage](#) (uint32_t NrChannel)
gets the final discharge voltage; unit: mV
- [uint32_t GetDischargeCapacity](#) (uint32_t NrChannel)
gets the discharge capacity; unit: μ Ah
- [uint32_t GetChargeCapacity](#) (uint32_t NrChannel)
gets the charge capacity; unit: μ Ah
- [uint32_t GetBatteryVoltage](#) (uint32_t NrChannel)
gets the battery voltage; unit: mV
- [uint32_t GetChannels](#) ()
gets number of channels
- void [SetRatedCapacityVolatile](#) (uint32_t NrChannel, [MbcRatedCapacityEnumNet](#) NewRatedCapacity)
sets the rated capacity (i.e. charge current) without storing it persistently
- void [SetChargingMode](#) (uint32_t NrChannel, [MbcChargingModeEnumNet](#) NewOperatingMode)

- sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge*
- [MbcChargingModeEnumNet GetChargingMode](#) (uint32_t NrChannel)
gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge
- [MbcChannelStateEnumNet GetChannelState](#) (uint32_t NrChannel)
gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge
- void [CapacityTest](#) (uint32_t NrChannel)
start capacity test on channel
- void [ChannelReset](#) (uint32_t NrChannel)
cancel charging and capacity test functions; check if battery is connected
- void [SetChargingPCoefficient](#) (uint32_t pCoefficient)
sets the p-coefficient for charging in mA/V / nominal charging current
- uint32_t [GetChargingPCoefficient](#) ()
gets the p-coefficient for charging in mA/V / nominal charging current
- void [SetRatedCapacity](#) (uint32_t NrChannel, [MbcRatedCapacityEnumNet](#) NewRatedCapacity)
sets the rated capacity
- [MbcRatedCapacityEnumNet GetRatedCapacity](#) (uint32_t NrChannel)
gets the rated capacity

Additional Inherited Members

11.73.1 Detailed Description

[CMultiBatteryChargerDeviceNet](#) is the class to access the MBC-08 device.

11.73.2 Constructor & Destructor Documentation

11.73.2.1 CMultiBatteryChargerDeviceNet() [CMultiBatteryChargerDeviceNet](#) ()

Initializes a new instance of the [CMultiBatteryChargerDeviceNet](#) class.

11.73.2.2 ~CMultiBatteryChargerDeviceNet() [virtual ~CMultiBatteryChargerDeviceNet](#) () [virtual]

11.73.2.3 "!CMultiBatteryChargerDeviceNet() [!CMultiBatteryChargerDeviceNet](#) ()

11.73.3 Member Function Documentation

11.73.3.1 CapacityTest() [void CapacityTest](#) (uint32_t NrChannel)

start capacity test on channel

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

11.73.3.2 ChannelReset() `void ChannelReset (`
 `uint32_t NrChannel)`

cancel charging and capacity test functions; check if battery is connected

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

11.73.3.3 GetBatteryVoltage() `uint32_t GetBatteryVoltage (`
 `uint32_t NrChannel)`

gets the battery voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the battery voltage in mV

11.73.3.4 GetChannels() `uint32_t GetChannels ()`

gets number of channels

Returns

number of channels

11.73.3.5 GetChannelState() `MbcChannelStateEnumNet GetChannelState (`
 `uint32_t NrChannel)`

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the current state

11.73.3.6 GetChargeCapacity() `uint32_t GetChargeCapacity (uint32_t NrChannel)`

gets the charge capacity; unit: μAh

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity in μAh

11.73.3.7 GetChargeCurrent() `uint32_t GetChargeCurrent (uint32_t NrChannel)`

gets the charge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the measured charge current in mA

11.73.3.8 GetChargingMode() `MbcChargingModeEnumNet GetChargingMode (uint32_t NrChannel)`

gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the charging mode

11.73.3.9 GetChargingPCoefficient() `uint32_t GetChargingPCoefficient ()`

gets the p-coefficient for charging in mA/V / nominal charging current

Returns

the p-coefficient

11.73.3.10 GetDischargeCapacity() `uint32_t GetDischargeCapacity (
uint32_t NrChannel)`

gets the discharge capacity; unit: μ Ah

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity in uAh

11.73.3.11 GetDischargeCurrent() `uint32_t GetDischargeCurrent (
uint32_t NrChannel)`

gets the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the measured discharge current in mA

11.73.3.12 GetDischargeCurrentSetPoint() `uint32_t GetDischargeCurrentSetPoint (
uint32_t NrChannel)`

gets the setpoint for the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the discharge current in mA

11.73.3.13 GetFinalDischargeVoltage() `uint32_t GetFinalDischargeVoltage (uint32_t NrChannel)`

gets the final discharge voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the battery voltage in mV at the end of discharge

11.73.3.14 GetRatedCapacity() `MbcRatedCapacityEnumNet GetRatedCapacity (uint32_t NrChannel)`

gets the rated capacity

Parameters

<i>NrChannel</i>	the channel number
------------------	--------------------

Returns

the capacity

11.73.3.15 SetChargingMode() `void SetChargingMode (uint32_t NrChannel, MbcChargingModeEnumNet NewOperatingMode)`

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

<i>NrChannel</i>	the channel number
<i>NewOperatingMode</i>	the charging mode

11.73.3.16 SetChargingPCoefficient() `void SetChargingPCoefficient (`
 `uint32_t pCoefficient)`

sets the p-coefficient for charging in mA/V / nominal charging current

Parameters

<i>pCoefficient</i>	the p-coefficient
---------------------	-------------------

11.73.3.17 SetDischargeCurrentSetPoint() `void SetDischargeCurrentSetPoint (`
 `uint32_t NrChannel,`
 `uint32_t DischargeCurrent_mA)`

sets the setpoint for the discharge current; unit: mA

Parameters

<i>NrChannel</i>	the channel number
<i>DischargeCurrent_mA</i>	the discharge current in mA

11.73.3.18 SetFinalDischargeVoltage() `void SetFinalDischargeVoltage (`
 `uint32_t NrChannel,`
 `uint32_t FinalDischargeVoltage_mV)`

sets the final discharge voltage; unit: mV

Parameters

<i>NrChannel</i>	the channel number
<i>FinalDischargeVoltage_mV</i>	the battery voltage in mV at the end of discharge

11.73.3.19 SetRatedCapacity() `void SetRatedCapacity (`
 `uint32_t NrChannel,`
 `MbcRatedCapacityEnumNet NewRatedCapacity)`

sets the rated capacity

Parameters

<i>NrChannel</i>	the channel number
<i>NewRatedCapacity</i>	the capacity

11.73.3.20 SetRatedCapacityVolatile() `void SetRatedCapacityVolatile (`
`uint32_t NrChannel,`
`MbcRatedCapacityEnumNet NewRatedCapacity)`

sets the rated capacity (i.e. charge current) without storing it persistently


Parameters

<i>NrChannel</i>	the channel number
<i>NewRatedCapacity</i>	the capacity

11.74 CMultiwellCallbackFunctionNet Class Reference

[CMultiwellCallbackFunctionNet](#) is the class to access the Multiwell-Mini-Stimulator

Inheritance diagram for [CMultiwellCallbackFunctionNet](#):

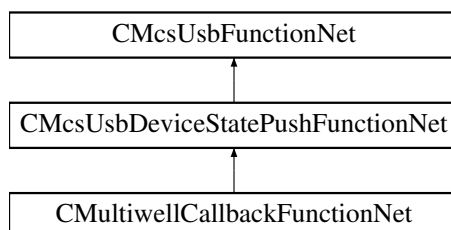

Public Member Functions

- delegate void [OnGetPlateClampStateByHeadstage](#) (uint32_t Headstage, [PlateClampEnumNet](#) plateState)
- [CMultiwellCallbackFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pMultiwellCallbackFunctionPointerContainer)
Initializes a new instance of the [CMultiwellCallbackFunctionNet](#) class.
- [CMultiwellCallbackFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual [~CMultiwellCallbackFunctionNet](#) ()
- [!CMultiwellCallbackFunctionNet](#) ()
- [PlateClampEnumNet](#) [GetPlateClampStateByHeadstage](#) (uint32_t Headstage)
Gets the state of the plate

Events

- [OnGetPlateClampStateByHeadstage](#)^ [GetPlateClampStateByHeadstageEvent](#) [add, remove, raise]
Event fires when the plate state for the headstage number has changed

Additional Inherited Members

11.74.1 Detailed Description

[CMultiwellCallbackFunctionNet](#) is the class to access the Multiwell-Mini-Stimulator

11.74.2 Constructor & Destructor Documentation

11.74.2.1 CMultiwellCallbackFunctionNet() [1/2] [CMultiwellCallbackFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *pMultiwellCallbackFunctionPointerContainer*)

Initializes a new instance of the [CMultiwellCallbackFunctionNet](#) class.

11.74.2.2 CMultiwellCallbackFunctionNet() [2/2] [CMultiwellCallbackFunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

11.74.2.3 ~CMultiwellCallbackFunctionNet() virtual [~CMultiwellCallbackFunctionNet](#) () [virtual]

11.74.2.4 !CMultiwellCallbackFunctionNet() [!CMultiwellCallbackFunctionNet](#) ()

11.74.3 Member Function Documentation

11.74.3.1 GetPlateClampStateByHeadstage() [PlateClampEnumNet](#) GetPlateClampStateByHeadstage (
 uint32_t *Headstage*)

Gets the state of the plate

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The plate state

11.74.3.2 OnGetPlateClampStateByHeadstage() `delegate void OnGetPlateClampStateByHeadstage (`
`uint32_t Headstage,`
`PlateClampEnumNet plateState)`

11.74.4 Event Documentation


11.74.4.1 GetPlateClampStateByHeadstageEvent `OnGetPlateClampStateByHeadstage^ GetPlateClamp↔`
`StateByHeadstageEvent [add], [remove], [raise]`

Event fires when the plate state for the headstage number has changed

11.75 CMultiwellDeviceNet Class Reference

[CMultiwellDeviceNet](#) is the class to access the Multiwell device.

Inheritance diagram for CMultiwellDeviceNet:

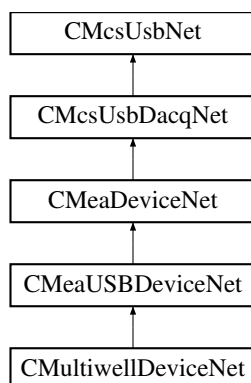

Public Member Functions

- [CMultiwellDeviceNet](#) ()
Initializes a new instance of the [CMultiwellDeviceNet](#) class.
- virtual [~CMultiwellDeviceNet](#) ()
- [!CMultiwellDeviceNet](#) ()
- [PlateClampEnumNet GetPlateClampState](#) ()
Gets the state of the Multiwell plate clamp.
- [PlateClampEnumNet GetPlateClampState](#) (uint32_t Headstage)
Gets the state of the plate
- void [OpenPlateClamp](#) ()
Opens the plate clamp.
- void [ClosePlateClamp](#) ()
Closes the plate clamp.

- void [StopPlateClamp](#) ()
Stops the plate clamp movement.
- uint32_t [GetPlateClampLockState](#) ()
Gets the state of the plate clamp lock.
- void [LockPlateClamp](#) ()
Locks the plate clamp.
- void [UnlockPlateClamp](#) ()
Unlocks the plate clamp.
- [MultiwellPlateTypeEnumNet](#) [GetPlateType](#) ()
Gets the plate type.
- [MultiwellPlateTypeEnumNet](#) [GetPlateType](#) (uint32_t Headstage)
Gets the plate type.
- void [SetPlateType](#) ([MultiwellPlateTypeEnumNet](#) plateType)
Sets the plate type.
- void [SetPlateType](#) (uint32_t Headstage, [MultiwellPlateTypeEnumNet](#) plateType)
Sets the plate type.
- void [SetPlateMux](#) (uint32_t muxSelection)
Selects a one quarter of the electrodes on a high density Multiwell plate.
- void [SetPlateMux](#) (uint32_t Headstage, uint32_t muxSelection)
Selects a one quarter of the electrodes on a high density Multiwell plate.
- uint32_t [GetPlateMux](#) ()
Gets the selected quarter of the electrodes on a high density Multiwell plate.
- uint32_t [GetPlateMux](#) (uint32_t Headstage)
Gets the selected quarter of the electrodes on a high density Multiwell plate.
- bool [IsPlateTypeValid](#) ()
Checks whether the plate type is valid, meaning all pins have contact.
- bool [IsPlateTypeValid](#) (uint32_t Headstage)
Checks whether the plate type is valid, meaning all pins have contact.
- void [SetPowerMuxPlate](#) (uint32_t Headstage, bool powerOn)
On the Multiwell Mini device, turn Power to the MUX Plate On or Off.
- bool [GetPowerMuxPlate](#) (uint32_t Headstage)
On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.
- void [SetTouchPadEnable](#) (uint32_t Headstage, bool Enable)
Enables or disables manual opening/closing of plate clamp via touch pad.
- bool [GetTouchPadEnable](#) (uint32_t Headstage)
Manual opening/closing of plate clamp via touch pad can be disabled.
- void [SetVolatileClampOffset](#) (uint32_t Headstage, int32_t CoverLipThickness_um)
The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.
- int32_t [GetVolatileClampOffset](#) (uint32_t Headstage)
The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Additional Inherited Members

11.75.1 Detailed Description

[CMultiwellDeviceNet](#) is the class to access the Multiwell device.

11.75.2 Constructor & Destructor Documentation

11.75.2.1 CMultiwellDeviceNet() `CMultiwellDeviceNet ()`

Initializes a new instance of the `CMultiwellDeviceNet` class.

11.75.2.2 ~CMultiwellDeviceNet() `virtual ~CMultiwellDeviceNet () [virtual]`**11.75.2.3 !CMultiwellDeviceNet()** `!CMultiwellDeviceNet ()`**11.75.3 Member Function Documentation****11.75.3.1 ClosePlateClamp()** `void ClosePlateClamp ()`

Closes the plate clamp.

11.75.3.2 GetPlateClampLockState() `uint32_t GetPlateClampLockState ()`

Gets the state of the plate clamp lock.

Returns

the state of the plate lock (unlocked/locked)

11.75.3.3 GetPlateClampState() [1/2] `PlateClampEnumNet GetPlateClampState ()`

Gets the state of the Multiwell plate clamp.

Returns

the state of the plate clamp (open/closed)

11.75.3.4 GetPlateClampState() [2/2] `PlateClampEnumNet GetPlateClampState (uint32_t Headstage)`

Gets the state of the plate

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The plate state

11.75.3.5 GetPlateMux() [1/2] `uint32_t GetPlateMux ()`

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Returns

the selected quarter

11.75.3.6 GetPlateMux() [2/2] `uint32_t GetPlateMux (
uint32_t Headstage)`

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

the selected quarter

11.75.3.7 GetPlateType() [1/2] `MultiwellPlateTypeEnumNet GetPlateType ()`

Gets the plate type.

Returns

the plate type

11.75.3.8 GetPlateType() [2/2] `MultiwellPlateTypeEnumNet GetPlateType (
uint32_t Headstage)`

Gets the plate type.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

the plate type

11.75.3.9 GetPowerMuxPlate() `bool GetPowerMuxPlate (`
`uint32_t Headstage)`

On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

"true" Power is On, "false" Power is Off

11.75.3.10 GetTouchPadEnable() `bool GetTouchPadEnable (`
`uint32_t Headstage)`

Manual opening/closing of plate clamp via touch pad can be disabled.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

"true" when plate clamp can be driven manually, otherwise "false".

11.75.3.11 GetVolatileClampOffset() `int32_t GetVolatileClampOffset (`
`uint32_t Headstage)`

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Parameters

<i>Headstage</i>	The headstage to be affected.
------------------	-------------------------------

Returns

Clamp distance reduction in um.

11.75.3.12 IsPlateTypeValid() [1/2] `bool IsPlateTypeValid ()`

Checks whether the plate type is valid, meaning all pins have contact.

Returns

"true" when all pins have contact, otherwise "false".

11.75.3.13 IsPlateTypeValid() [2/2] `bool IsPlateTypeValid (
 uint32_t Headstage)`

Checks whether the plate type is valid, meaning all pins have contact.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

"true" when all pins have contact, otherwise "false".

11.75.3.14 LockPlateClamp() `void LockPlateClamp ()`

Locks the plate clamp.

11.75.3.15 OpenPlateClamp() `void OpenPlateClamp ()`

Opens the plate clamp.

11.75.3.16 SetPlateMux() [1/2] `void SetPlateMux (
 uint32_t Headstage,
 uint32_t muxSelection)`

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>Headstage</i>	The headstage to query.
<i>muxSelection</i>	the selected quarter

11.75.3.17 SetPlateMux() [2/2] `void SetPlateMux (`
`uint32_t muxSelection)`

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

<i>muxSelection</i>	the selected quarter
---------------------	----------------------

11.75.3.18 SetPlateType() [1/2] `void SetPlateType (`
`MultiwellPlateTypeEnumNet plateType)`

Sets the plate type.

Parameters

<i>plateType</i>	the plate type
------------------	----------------

11.75.3.19 SetPlateType() [2/2] `void SetPlateType (`
`uint32_t Headstage,`
`MultiwellPlateTypeEnumNet plateType)`

Sets the plate type.

Parameters

<i>Headstage</i>	The headstage to query.
<i>plateType</i>	the plate type

11.75.3.20 SetPowerMuxPlate() `void SetPowerMuxPlate (`
`uint32_t Headstage,`
`bool powerOn)`

On the Multiwell Mini device, turn Power to the MUX Plate On or Off.

Parameters

<i>Headstage</i>	The headstage to query.
<i>powerOn</i>	"true" to turn Power On, "false" to turn Power Off

11.75.3.21 SetTouchPadEnable() `void SetTouchPadEnable (`
 `uint32_t Headstage,`
 `bool Enable)`

Enables or disables manual opening/closing of plate clamp via touch pad.

Parameters

<i>Headstage</i>	The headstage to be affected.
<i>Enable</i>	"true" when plate clamp shall be driven manually, otherwise "false".

11.75.3.22 SetVolatileClampOffset() `void SetVolatileClampOffset (`
 `uint32_t Headstage,`
 `int32_t CoverLipThickness_um)`

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Parameters

<i>Headstage</i>	The headstage to be affected.
<i>CoverLipThickness_um</i>	Clamp distance reduction in um. Range: 200um..400um; Typical value: 300um

11.75.3.23 StopPlateClamp() `void StopPlateClamp ()`

Stops the plate clamp movement.


11.75.3.24 UnlockPlateClamp() `void UnlockPlateClamp ()`

Unlocks the plate clamp.

11.76 CMultiwellOptoStimFunctionNet Class Reference

[CMultiwellOptoStimFunctionNet](#) is the class to access the optical properties of the Multiwell Optostim device

Inheritance diagram for [CMultiwellOptoStimFunctionNet](#):

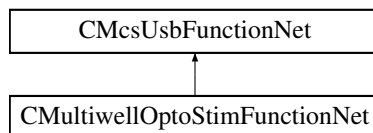

Public Member Functions

- [CMultiwellOptoStimFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pMultiwellOptoStimFunctionPointerContainer)
Initializes a new instance of the [CMultiwellOptoStimFunctionNet](#) class.
- [CMultiwellOptoStimFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CMultiwellOptoStimFunctionNet](#) ()
- [!CMultiwellOptoStimFunctionNet](#) ()
- [uint32_t](#) [GetWaveLengthInNanometer](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetAbsMaxCurrentInMicroAmp](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetMaxDurationHighCurrentInMicroSec](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetMaxDutyCycleHighCurrent](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetPermanentCurrentInMicroAmp](#) ([uint16_t](#) channel)
- [uint32_t](#) [GetColorRgb](#) ([uint16_t](#) channel)
- [String](#)[^] [GetColorStr](#) ([uint16_t](#) channel)
- void [SetWaveLengthInNanometer](#) ([uint16_t](#) channel, [uint32_t](#) WaveLength_nm)
- void [SetAbsMaxCurrentInMicroAmp](#) ([uint16_t](#) channel, [uint32_t](#) AbsoluteMaxCurrent_uA)
- void [SetMaxDurationHighCurrentInMicroSec](#) ([uint16_t](#) channel, [uint32_t](#) AbsoluteMaxDuration_us)
- void [SetMaxDutyCycleHighCurrent](#) ([uint16_t](#) channel, [uint32_t](#) MaxDutyCycleHighCurrent)
- void [SetPermanentCurrentInMicroAmp](#) ([uint16_t](#) channel, [uint32_t](#) PermanentCurrent_uA)
- void [SetColorRgb](#) ([uint16_t](#) channel, [uint32_t](#) ColorRGB)
- void [SetColorStr](#) ([uint16_t](#) channel, [String](#)[^] ColorString)

Additional Inherited Members

11.76.1 Detailed Description

[CMultiwellOptoStimFunctionNet](#) is the class to access the optical properties of the Multiwell Optostim device

11.76.2 Constructor & Destructor Documentation

11.76.2.1 [CMultiwellOptoStimFunctionNet](#)() [1/2] [CMultiwellOptoStimFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pMultiwellOptoStimFunctionPointerContainer)

Initializes a new instance of the [CMultiwellOptoStimFunctionNet](#) class.

11.76.2.2 CMultiwellOptoStimFunctionNet() [2/2] `CMultiwellOptoStimFunctionNet (CMcsUsbNet^ mcsusb)`

11.76.2.3 ~CMultiwellOptoStimFunctionNet() `virtual ~CMultiwellOptoStimFunctionNet ()` [virtual]

11.76.2.4 ~!CMultiwellOptoStimFunctionNet() `!CMultiwellOptoStimFunctionNet ()`

11.76.3 Member Function Documentation

11.76.3.1 GetAbsMaxCurrentInMicroAmp() `uint32_t GetAbsMaxCurrentInMicroAmp (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

absolute max. current; unit: uA

11.76.3.2 GetColorRgb() `uint32_t GetColorRgb (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

RGB-value of LED color

11.76.3.3 GetColorStr() `String ^ GetColorStr (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

LED color as string

11.76.3.4 GetMaxDurationHighCurrentInMicroSec() `uint32_t GetMaxDurationHighCurrentInMicroSec (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. duration the LED can stand the abs. max current; unit: us

11.76.3.5 GetMaxDutyCycleHighCurrent() `uint32_t GetMaxDutyCycleHighCurrent (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. duty cycle at max. current; unit: 100*%

11.76.3.6 GetPermanentCurrentInMicroAmp() `uint32_t GetPermanentCurrentInMicroAmp (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

max. current the LED can stand when always switched on; unit: uA

11.76.3.7 GetWaveLengthInNanometer() `uint32_t GetWaveLengthInNanometer (uint16_t channel)`

Parameters

<i>channel</i>	the (analog) channel number
----------------	-----------------------------

Returns

wavelength of this channel's LEDs; unit: nm

11.76.3.8 SetAbsMaxCurrentInMicroAmp() `void SetAbsMaxCurrentInMicroAmp (uint16_t channel, uint32_t AbsoluteMaxCurrent_uA)`

Parameters

<i>channel</i>	the (analog) channel number
<i>AbsoluteMaxCurrent_uA</i>	absolute max. current; unit: uA

11.76.3.9 SetColorRgb() `void SetColorRgb (uint16_t channel, uint32_t ColorRGB)`

Parameters

<i>channel</i>	the (analog) channel number
<i>ColorRGB</i>	RGB-value of LED color

11.76.3.10 SetColorStr() `void SetColorStr (uint16_t channel, String^ ColorString)`

Parameters

<i>channel</i>	the (analog) channel number
<i>ColorString</i>	LED color as string

11.76.3.11 SetMaxDurationHighCurrentInMicroSec() `void SetMaxDurationHighCurrentInMicroSec (`

```
uint16_t channel,
uint32_t AbsoluteMaxDuration_us )
```

Parameters

<i>channel</i>	the (analog) channel number
<i>AbsoluteMaxDuration_us</i>	max. duration the LED can stand the abs. max current; unit: us

11.76.3.12 SetMaxDutyCycleHighCurrent() `void SetMaxDutyCycleHighCurrent (`

```
uint16_t channel,
uint32_t MaxDutyCycleHighCurrent )
```

Parameters

<i>channel</i>	the (analog) channel number
<i>MaxDutyCycleHighCurrent</i>	max. duty cycle at max. current; unit: 100*%

11.76.3.13 SetPermanentCurrentInMicroAmp() `void SetPermanentCurrentInMicroAmp (`

```
uint16_t channel,
uint32_t PermanentCurrent_uA )
```

Parameters

<i>channel</i>	the (analog) channel number
<i>PermanentCurrent_uA</i>	max. current the LED can stand when always switched on; unit: uA

11.76.3.14 SetWaveLengthInNanometer() `void SetWaveLengthInNanometer (`


```
uint16_t channel,
uint32_t WaveLength_nm )
```

Parameters

<i>channel</i>	the (analog) channel number
<i>WaveLength_nm</i>	wavelength of this channel's LEDs; unit: nm

11.77 CNF_GenDeviceNet Class Reference

Inheritance diagram for CNF_GenDeviceNet:

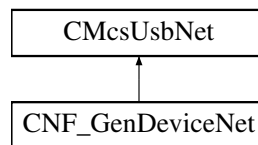

Public Member Functions

- [CNF_GenDeviceNet](#) (void)
- [~CNF_GenDeviceNet](#) (void)
- void [Set_Values](#) (unsigned int frequency, unsigned int amplitude)

Additional Inherited Members

11.77.1 Constructor & Destructor Documentation

11.77.1.1 [CNF_GenDeviceNet\(\)](#) [CNF_GenDeviceNet](#) (
 void)


11.77.1.2 [~CNF_GenDeviceNet\(\)](#) [~CNF_GenDeviceNet](#) (
 void)

11.77.2 Member Function Documentation

11.77.2.1 [Set_Values\(\)](#) void [Set_Values](#) (
 unsigned int *frequency*,
 unsigned int *amplitude*)

11.78 COctoPotDeviceNet Class Reference

Inheritance diagram for COctoPotDeviceNet:

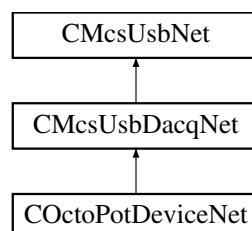

Public Member Functions

- [COctoPotDeviceNet](#) (void)
- [COctoPotDeviceNet](#) ([OnChannelData](#)^ channelData, [OnError](#)^ error)
- [uint32_t SetOutputRate](#) ([uint32_t](#) rate)
- [uint32_t SetBathclamp](#) (unsigned int block, bool enable)
- [uint32_t SetDacValue](#) (int channel, int value)
- [uint32_t SetDacAutoControl](#) (unsigned int channel)
- [uint32_t SetPidParameter](#) (unsigned int channel, int const_p, int const_i, int shift_p, int shift_i)
- [uint32_t SetRampParameter](#) (unsigned int channel, int start, int min, int max, int slope, int slope2, int pause, unsigned int samples)
- [uint32_t RampStart](#) (int channelmap)
- [uint32_t SetSineParameter](#) (unsigned int channel, int amplitude)
- [uint32_t SineStart](#) (int channelmap)
- [uint32_t SetPatternListEntry](#) (unsigned int channel, unsigned int position, unsigned int duration, int value)
- [uint32_t PatternListStart](#) (int channelmap)
- [uint32_t SetAdcOffset](#) (unsigned int channel, int offset)
- [uint32_t SetDacOffset](#) (unsigned int channel, int offset)
- [uint32_t ResetAdcOffset](#) (unsigned int channel)
- [uint32_t ResetDacOffset](#) (unsigned int channel)
- [uint32_t BurnAdcOffset](#) ()
- [uint32_t BurnDacOffset](#) ()
- [uint32_t GetAdcOffset](#) (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- [uint32_t GetDacOffset](#) (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- [uint32_t SetAmplificationSwitch](#) (unsigned int channel, unsigned int state)
- [uint32_t SetChannelSwitch](#) (unsigned int channel, unsigned int state)
- [uint32_t SetNumberOfChannels](#) (unsigned int NumberOfChannels)
- [uint32_t EnableDigitalIn](#) (bool enable)
- [uint32_t EnableTimestamp](#) (bool enable)
- [uint32_t EnableChecksum](#) (bool enable)

Additional Inherited Members

11.78.1 Constructor & Destructor Documentation

11.78.1.1 COctoPotDeviceNet() [1/2] [COctoPotDeviceNet](#) (
 void)

11.78.1.2 COctoPotDeviceNet() [2/2] [COctoPotDeviceNet](#) (
 [OnChannelData](#)^ channelData,
 [OnError](#)^ error)

11.78.2 Member Function Documentation

11.78.2.1 BurnAdcOffset() uint32_t BurnAdcOffset ()

11.78.2.2 BurnDacOffset() uint32_t BurnDacOffset ()

11.78.2.3 EnableChecksum() uint32_t EnableChecksum (
bool *enable*)

11.78.2.4 EnableDigitalIn() uint32_t EnableDigitalIn (
bool *enable*)

11.78.2.5 EnableTimestamp() uint32_t EnableTimestamp (
bool *enable*)

11.78.2.6 GetAdcOffset() uint32_t GetAdcOffset (
unsigned int *channel*,
[System::Runtime::InteropServices::Out] int ^ *offset*)

11.78.2.7 GetDacOffset() uint32_t GetDacOffset (
unsigned int *channel*,
[System::Runtime::InteropServices::Out] int ^ *offset*)

11.78.2.8 PatternListStart() uint32_t PatternListStart (
int *channelmap*)

11.78.2.9 RampStart() uint32_t RampStart (
int *channelmap*)

11.78.2.10 ResetAdcOffset() uint32_t ResetAdcOffset (
unsigned int *channel*)

11.78.2.11 ResetDacOffset() uint32_t ResetDacOffset (
 unsigned int *channel*)

11.78.2.12 SetAdcOffset() uint32_t SetAdcOffset (
 unsigned int *channel*,
 int *offset*)

11.78.2.13 SetAmplificationSwitch() uint32_t SetAmplificationSwitch (
 unsigned int *channel*,
 unsigned int *state*)

11.78.2.14 SetBathclamp() uint32_t SetBathclamp (
 unsigned int *block*,
 bool *enable*)

11.78.2.15 SetChannelSwitch() uint32_t SetChannelSwitch (
 unsigned int *channel*,
 unsigned int *state*)

11.78.2.16 SetDacAutoControl() uint32_t SetDacAutoControl (
 unsigned int *channel*)

11.78.2.17 SetDacOffset() uint32_t SetDacOffset (
 unsigned int *channel*,
 int *offset*)

11.78.2.18 SetDacValue() uint32_t SetDacValue (
 int *channel*,
 int *value*)

11.78.2.19 SetNumberOfChannels() uint32_t SetNumberOfChannels (
 unsigned int *NumberOfChannels*)

11.78.2.20 SetOutputRate() uint32_t SetOutputRate (
uint32_t *rate*)

11.78.2.21 SetPatternListEntry() uint32_t SetPatternListEntry (
unsigned int *channel*,
unsigned int *position*,
unsigned int *duration*,
int *value*)

11.78.2.22 SetPidParameter() uint32_t SetPidParameter (
unsigned int *channel*,
int *const_p*,
int *const_i*,
int *shift_p*,
int *shift_i*)


11.78.2.23 SetRampParameter() uint32_t SetRampParameter (
unsigned int *channel*,
int *start*,
int *min*,
int *max*,
int *slope*,
int *slope2*,
int *pause*,
unsigned int *samples*)

11.78.2.24 SetSineParameter() uint32_t SetSineParameter (
unsigned int *channel*,
int *amplitude*)

11.78.2.25 SineStart() uint32_t SineStart (
int *channelmap*)

11.79 COkuvisionStimulatorDeviceNet Class Reference

Inheritance diagram for COkuvisionStimulatorDeviceNet:

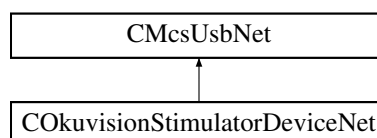

Public Member Functions

- [COKuvisionStimulatorDeviceNet](#) (void)
- [~COKuvisionStimulatorDeviceNet](#) (void)
- void [SetPulseform](#) (int channel, int current, int pulsewidth, int periode, int duration)
- void [GetPulseform](#) (int channel, [System::Runtime::InteropServices::Out] int% current, [System::Runtime::InteropServices::Out] int% pulsewidth, [System::Runtime::InteropServices::Out] int% periode, [System::Runtime::InteropServices::Out] int% duration)
- void [SetMaxPower](#) (int channel, int power)
- int [GetMaxPower](#) (int channel)
- void [SetMaxVoltage](#) (int channel, int voltage)
- int [GetMaxVoltage](#) (int channel)
- void [SetCheckVoltage](#) (int channel, int voltage)
- int [GetCheckVoltage](#) (int channel)
- int [GetVoltage](#) (int channel)
- void [SetDACOffset](#) (int channel, int part, int offset)
- int [GetDACOffset](#) (int channel, int part)
- void [SetRTC](#) (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)
- void [GetRTC](#) ([System::Runtime::InteropServices::Out] uint8_t% year, [System::Runtime::InteropServices::Out] uint8_t% month, [System::Runtime::InteropServices::Out] uint8_t% day, [System::Runtime::InteropServices::Out] uint8_t% hour, [System::Runtime::InteropServices::Out] uint8_t% minute, [System::Runtime::InteropServices::Out] uint8_t% second)
- void [SetRTC](#) (DateTime timestamp)
- DateTime [GetRTC](#) ()
- void [GetStimulatorStatus](#) ([System::Runtime::InteropServices::Out] int% startstop, [System::Runtime::InteropServices::Out] int% last_error, [System::Runtime::InteropServices::Out] int% battery_status)
- void [SetCurrentFactor](#) (int channel, int factor)
- int [GetCurrentFactor](#) (int channel)

Additional Inherited Members

11.79.1 Constructor & Destructor Documentation

11.79.1.1 COKuvisionStimulatorDeviceNet() [COKuvisionStimulatorDeviceNet](#) (void)

11.79.1.2 ~COKuvisionStimulatorDeviceNet() [~COKuvisionStimulatorDeviceNet](#) (void)

11.79.2 Member Function Documentation

11.79.2.1 GetCheckVoltage() int GetCheckVoltage (int channel)

11.79.2.2 GetCurrentFactor() `int GetCurrentFactor (`
`int channel)`

11.79.2.3 GetDACOffset() `int GetDACOffset (`
`int channel,`
`int part)`

11.79.2.4 GetMaxPower() `int GetMaxPower (`
`int channel)`

11.79.2.5 GetMaxVoltage() `int GetMaxVoltage (`
`int channel)`

11.79.2.6 GetPulseform() `void GetPulseform (`
`int channel,`
`[System::Runtime::InteropServices::Out] int% current,`
`[System::Runtime::InteropServices::Out] int% pulsewidth,`
`[System::Runtime::InteropServices::Out] int% periode,`
`[System::Runtime::InteropServices::Out] int% duration)`

11.79.2.7 GetRTC() `[1/2] DateTime GetRTC ()`

11.79.2.8 GetRTC() `[2/2] void GetRTC (`
`[System::Runtime::InteropServices::Out] uint8_t% year,`
`[System::Runtime::InteropServices::Out] uint8_t% month,`
`[System::Runtime::InteropServices::Out] uint8_t% day,`
`[System::Runtime::InteropServices::Out] uint8_t% hour,`
`[System::Runtime::InteropServices::Out] uint8_t% minute,`
`[System::Runtime::InteropServices::Out] uint8_t% second)`

11.79.2.9 GetStimulatorStatus() `void GetStimulatorStatus (`
`[System::Runtime::InteropServices::Out] int% startstop,`
`[System::Runtime::InteropServices::Out] int% last_error,`
`[System::Runtime::InteropServices::Out] int% battery_status)`

11.79.2.10 GetVoltage() `int GetVoltage (`
`int channel)`

11.79.2.11 SetCheckVoltage() `void SetCheckVoltage (`
`int channel,`
`int voltage)`

11.79.2.12 SetCurrentFactor() `void SetCurrentFactor (`
`int channel,`
`int factor)`

11.79.2.13 SetDACOffset() `void SetDACOffset (`
`int channel,`
`int part,`
`int offset)`

11.79.2.14 SetMaxPower() `void SetMaxPower (`
`int channel,`
`int power)`

11.79.2.15 SetMaxVoltage() `void SetMaxVoltage (`
`int channel,`
`int voltage)`

11.79.2.16 SetPulseform() `void SetPulseform (`
`int channel,`
`int current,`
`int pulsewidth,`
`int periode,`
`int duration)`


11.79.2.17 SetRTC() `[1/2] void SetRTC (`
`DateTime timestamp)`

11.79.2.18 SetRTC() [2/2] `void SetRTC (`
`uint8_t year,`
`uint8_t month,`
`uint8_t day,`
`uint8_t hour,`
`uint8_t minute,`
`uint8_t second)`

11.80 CPatchServerDeviceNet Class Reference

[CPatchServerDeviceNet](#) is the class to control the MCS PatchServer device

Inheritance diagram for CPatchServerDeviceNet:

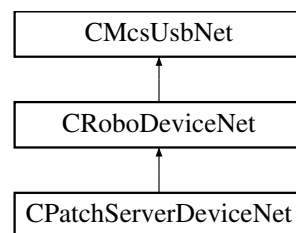

Public Member Functions

- [CPatchServerDeviceNet](#) (void)

Properties

- [CMcsBus_SensorNet](#)[^] [Sensor](#) [get]

Additional Inherited Members

11.80.1 Detailed Description

[CPatchServerDeviceNet](#) is the class to control the MCS PatchServer device

11.80.2 Constructor & Destructor Documentation


11.80.2.1 CPatchServerDeviceNet() [CPatchServerDeviceNet](#) (
`void`)

11.80.3 Property Documentation

11.80.3.1 Sensor `CMcsBus_SensorNet^ Sensor [get]`

11.81 CPathIdentDeviceNet Class Reference

Inheritance diagram for CPathIdentDeviceNet:

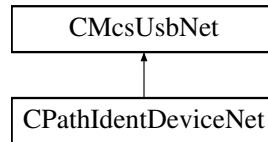

Public Member Functions

- `CPathIdentDeviceNet` (void)
- `~CPathIdentDeviceNet` (void)
- void `Set_Values` (unsigned int frequency, unsigned int amplitude)
- void `Measure` ([System::Runtime::InteropServices::Out] unsigned int% phase, [System::Runtime::InteropServices::Out] unsigned int% amplitude)

Additional Inherited Members

11.81.1 Constructor & Destructor Documentation

11.81.1.1 CPathIdentDeviceNet() `CPathIdentDeviceNet (void)`

11.81.1.2 ~CPathIdentDeviceNet() `~CPathIdentDeviceNet (void)`


11.81.2 Member Function Documentation

11.81.2.1 Measure() `void Measure ([System::Runtime::InteropServices::Out] unsigned int% phase, [System::Runtime::InteropServices::Out] unsigned int% amplitude)`

11.81.2.2 Set_Values() void Set_Values (
 unsigned int *frequency*,
 unsigned int *amplitude*)

11.82 CPedoterDeviceNet Class Reference

Inheritance diagram for CPedoterDeviceNet:

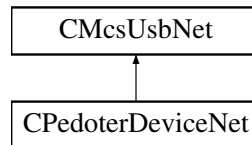

Public Member Functions

- [CPedoterDeviceNet](#) ()
Initializes a new instance of the [CPedoterDeviceNet](#) class.
- virtual [~CPedoterDeviceNet](#) ()
- [!CPedoterDeviceNet](#) ()
- uint32_t [GetCommand](#) (uint16_t Argument)
Get value from the pedoter device
- void [SetCommand](#) (uint16_t Argument, uint32_t pData)
Set value on the pedoter device

Additional Inherited Members

11.82.1 Detailed Description

11.82.2 Constructor & Destructor Documentation

11.82.2.1 CPedoterDeviceNet() [CPedoterDeviceNet](#) ()

Initializes a new instance of the [CPedoterDeviceNet](#) class.

11.82.2.2 ~CPedoterDeviceNet() virtual [~CPedoterDeviceNet](#) () [virtual]

11.82.2.3 "!CPedoterDeviceNet() [!CPedoterDeviceNet](#) ()

11.82.3 Member Function Documentation

11.82.3.1 GetCommand() uint32_t GetCommand (
 uint16_t *Argument*)

Get value from the pedoter device

Parameters

<i>Argument</i>	argument
-----------------	----------

Returns

value

11.82.3.2 SetCommand() void SetCommand (
 uint16_t *Argument*,
 uint32_t *pData*)

Set value on the pedoter device


Parameters

<i>Argument</i>	argument
<i>pData</i>	value

11.83 CPeristalticPumpDeviceNet Class Reference

[CPeristalticPumpDeviceNet](#) is the class to control a Persistaltic Pump.

Inheritance diagram for CPeristalticPumpDeviceNet:

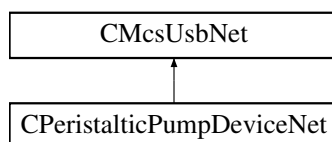**Public Member Functions**

- [CPeristalticPumpDeviceNet](#) (void)
 Initialize a new instance of the [CPeristalticPumpDeviceNet](#) class.
- [~CPeristalticPumpDeviceNet](#) (void)

Properties

- [CMcsBus_MotorControlNet](#)^ [McsBus_MotorControl](#) [get]

Additional Inherited Members**11.83.1 Detailed Description**

[CPeristalticPumpDeviceNet](#) is the class to control a Persistaltic Pump.

11.83.2 Constructor & Destructor Documentation

11.83.2.1 CPeristalticPumpDeviceNet() `CPeristalticPumpDeviceNet (void)`

Initialize a new instance of the `CPeristalticPumpDeviceNet` class.


11.83.2.2 ~CPeristalticPumpDeviceNet() `~CPeristalticPumpDeviceNet (void)`

11.83.3 Property Documentation

11.83.3.1 McsBus_MotorControl `CMcsBus_MotorControlNet^ McsBus_MotorControl [get]`

11.84 CPgaDeviceNet Class Reference

Inheritance diagram for CPgaDeviceNet:

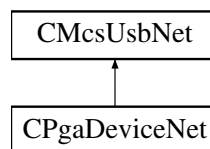

Public Member Functions

- `CPgaDeviceNet ()`
- `~CPgaDeviceNet ()`
- `uint32_t GetNumFrequencyRanges ([System::Runtime::InteropServices::Out]int% numRanges)`
- `uint32_t GetFrequencyRange (int rangeIndex, [System::Runtime::InteropServices::Out]int% low, [System::Runtime::InteropServices::Out]int% high, [System::Runtime::InteropServices::Out]int% channels, [System::Runtime::InteropServices::Out]int% gain)`
- `uint32_t GetNumAmplifications ([System::Runtime::InteropServices::Out]int% number)`
- `uint32_t GetAmplification (int index, [System::Runtime::InteropServices::Out]int% amplification, [System::Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)`
- `uint32_t DefineNumFrequencyRanges (int rnum)`
- `uint32_t DefineFrequencyRange (int index, int low, int high, int channels, int gain)`
- `uint32_t DefineNumAmplifications (int number)`
- `uint32_t DefineAmplification (int index, int amplification, int poti1, int poti2)`
- `uint32_t SetGain (int channel, int Gain, int poti1, int poti2)`
- `uint32_t GetGain (int channel, [System::Runtime::InteropServices::Out]int% Gain, [System::Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)`
- `uint32_t ApplyGains ()`

Additional Inherited Members

11.84.1 Constructor & Destructor Documentation

11.84.1.1 CPgaDeviceNet() `CPgaDeviceNet ()`

11.84.1.2 ~CPgaDeviceNet() `~CPgaDeviceNet ()`

11.84.2 Member Function Documentation

11.84.2.1 ApplyGains() `uint32_t ApplyGains ()`

11.84.2.2 DefineAmplification() `uint32_t DefineAmplification (`
 `int index,`
 `int amplification,`
 `int poti1,`
 `int poti2)`

11.84.2.3 DefineFrequencyRange() `uint32_t DefineFrequencyRange (`
 `int index,`
 `int low,`
 `int high,`
 `int channels,`
 `int gain)`

11.84.2.4 DefineNumAmplifications() `uint32_t DefineNumAmplifications (`
 `int number)`

11.84.2.5 DefineNumFrequencyRanges() `uint32_t DefineNumFrequencyRanges (`
 `int rnum)`

11.84.2.6 GetAmplification() uint32_t GetAmplification (
 int *index*,
 [System::Runtime::InteropServices::Out] int% *amplification*,
 [System::Runtime::InteropServices::Out] int% *poti1*,
 [System::Runtime::InteropServices::Out] int% *poti2*)

11.84.2.7 GetFrequencyRange() uint32_t GetFrequencyRange (
 int *rangeIndex*,
 [System::Runtime::InteropServices::Out] int% *low*,
 [System::Runtime::InteropServices::Out] int% *high*,
 [System::Runtime::InteropServices::Out] int% *channels*,
 [System::Runtime::InteropServices::Out] int% *gain*)

11.84.2.8 GetGain() uint32_t GetGain (
 int *channel*,
 [System::Runtime::InteropServices::Out] int% *Gain*,
 [System::Runtime::InteropServices::Out] int% *poti1*,
 [System::Runtime::InteropServices::Out] int% *poti2*)

11.84.2.9 GetNumAmplifications() uint32_t GetNumAmplifications (
 [System::Runtime::InteropServices::Out] int% *number*)


11.84.2.10 GetNumFrequencyRanges() uint32_t GetNumFrequencyRanges (
 [System::Runtime::InteropServices::Out] int% *numRanges*)

11.84.2.11 SetGain() uint32_t SetGain (
 int *channel*,
 int *Gain*,
 int *poti1*,
 int *poti2*)

11.85 CPositionIIDeviceNet Class Reference

[CPositionIIDeviceNet](#) is the class to control PositionII devices

Inheritance diagram for CPositionIIDeviceNet:

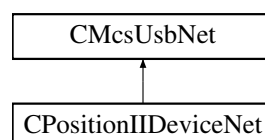

Public Member Functions

- [CPositionIIDeviceNet](#) ()
Initializes a new instance of the [CPositionIIDeviceNet](#) class.
- virtual [~CPositionIIDeviceNet](#) ()
- [!CPositionIIDeviceNet](#) ()
- uint32_t [GetCoilCommunication](#) (uint16_t coil)
get if the communication to the coil is working
- uint32_t [GetOnOff](#) (uint16_t coil)
get if the coil is switched on/off
- void [SwitchOnOff](#) (uint16_t coil, uint32_t on)
switched the coil on of
- uint32_t [GetImplantState](#) (uint16_t coil)
gets the implantat state
- uint32_t [GetImplantCurrentSetpoint](#) (uint16_t coil)
sets the implant current setpoint
- void [SetImplantCurrentSetpoint](#) (uint16_t coil, uint32_t current)
gets the implant current setpoint
- uint32_t [GetPowerStrength](#) (uint16_t coil)
sets the power for the trigger pulses
- void [SetPowerStrength](#) (uint16_t coil, uint32_t power)
gets the power for the trigger pulses
- uint32_t [GetImplantResult](#) (uint16_t coil)
gets the last result of the implant pulse trigger
- void [GetRTC](#) ([System::Runtime::InteropServices::Out]uint8_t% year, [System::Runtime::InteropServices::Out]uint8_t% month, [System::Runtime::InteropServices::Out]uint8_t% day, [System::Runtime::InteropServices::Out]uint8_t% hour, [System::Runtime::InteropServices::Out]uint8_t% minute, [System::Runtime::InteropServices::Out]uint8_t% second)
Get the RTC
- void [SetRTC](#) (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)
Set the RTC
- uint32_t [GetStateDebugData](#) (uint16_t coil)
get the debug queue state
- void [SetStateDebugData](#) (uint16_t coil, uint32_t state)
clears/starts/stops the debug queue for a certain coil
- void [GetDebugData](#) (uint16_t coil, [System::Runtime::InteropServices::Out]uint16_t% index, [System::Runtime::InteropServices::Out]uint16_t% voltage, [System::Runtime::InteropServices::Out]uint16_t% numberofpulses, [System::Runtime::InteropServices::Out]uint16_t% mediantime)
get the oldest debug entry for a certain coil
- uint32_t [GetStateEventData](#) ()
get the event queue state
- void [SetStateEventData](#) (uint32_t state)
clears/starts/stops the event queue for a certain coil
- void [GetEventData](#) ([System::Runtime::InteropServices::Out]uint16_t% index, [System::Runtime::InteropServices::Out]uint8_t% year, [System::Runtime::InteropServices::Out]uint8_t% month, [System::Runtime::InteropServices::Out]uint8_t% day, [System::Runtime::InteropServices::Out]uint8_t% hour, [System::Runtime::InteropServices::Out]uint8_t% minute, [System::Runtime::InteropServices::Out]uint8_t% second, [System::Runtime::InteropServices::Out]uint16_t% coil, [System::Runtime::InteropServices::Out]uint16_t% type, [System::Runtime::InteropServices::Out]uint16_t% value)
get the oldest event entry

Properties

- [CRFFunctionNet](#)[^] [RFFunction](#) [get]

Additional Inherited Members

11.85.1 Detailed Description

[CPositionIIDeviceNet](#) is the class to control PositionII devices

11.85.2 Constructor & Destructor Documentation

11.85.2.1 CPositionIIDeviceNet() [CPositionIIDeviceNet](#) ()

Initializes a new instance of the [CPositionIIDeviceNet](#) class.

11.85.2.2 ~CPositionIIDeviceNet() [virtual](#) [~CPositionIIDeviceNet](#) () [virtual]

11.85.2.3 "!CPositionIIDeviceNet() [!CPositionIIDeviceNet](#) ()

11.85.3 Member Function Documentation

11.85.3.1 GetCoilCommunication() [uint32_t](#) [GetCoilCommunication](#) ([uint16_t](#) *coil*)

get if the communication to the coil is working

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

is communicating

11.85.3.2 GetDebugData() void GetDebugData (

```

    uint16_t coil,
    [System::Runtime::InteropServices::Out] uint16_t% index,
    [System::Runtime::InteropServices::Out] uint16_t% voltage,
    [System::Runtime::InteropServices::Out] uint16_t% numberofpulses,
    [System::Runtime::InteropServices::Out] uint16_t% mediantime )

```

get the oldest debug entry for a certain coil

Parameters

<i>coil</i>	the coil
<i>index</i>	the debug entry index number
<i>voltage</i>	the voltage applied
<i>numberofpulses</i>	the number of pulses detected
<i>mediantime</i>	the median time between pulses

11.85.3.3 GetEventData() void GetEventData (

```

    [System::Runtime::InteropServices::Out] uint16_t% index,
    [System::Runtime::InteropServices::Out] uint8_t% year,
    [System::Runtime::InteropServices::Out] uint8_t% month,
    [System::Runtime::InteropServices::Out] uint8_t% day,
    [System::Runtime::InteropServices::Out] uint8_t% hour,
    [System::Runtime::InteropServices::Out] uint8_t% minute,
    [System::Runtime::InteropServices::Out] uint8_t% second,
    [System::Runtime::InteropServices::Out] uint16_t% coil,
    [System::Runtime::InteropServices::Out] uint16_t% type,
    [System::Runtime::InteropServices::Out] uint16_t% value )

```

get the oldest event entry

Parameters

<i>index</i>	the event index number
<i>year</i>	the year
<i>month</i>	the month
<i>day</i>	the day
<i>hour</i>	the hour
<i>minute</i>	the minute
<i>second</i>	the second
<i>coil</i>	the coil
<i>type</i>	the event type
<i>value</i>	the even value

11.85.3.4 GetImplantCurrentSetpoint() uint32_t GetImplantCurrentSetpoint (

```

    uint16_t coil )

```

sets the implant current setpoint

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the current

11.85.3.5 GetImplantResult() `uint32_t GetImplantResult (`
`uint16_t coil)`

gets the last result of the implant pulse trigger

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the result

11.85.3.6 GetImplantState() `uint32_t GetImplantState (`
`uint16_t coil)`

gets the implantat state

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the state

11.85.3.7 GetOnOff() `uint32_t GetOnOff (`
`uint16_t coil)`

get if the coil is switched on/off

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

0 = off, 1 = on

11.85.3.8 GetPowerStrength() `uint32_t GetPowerStrength (`
`uint16_t coil)`

sets the power for the trigger pulses

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the power

11.85.3.9 GetRTC() `void GetRTC (`
`[System::Runtime::InteropServices::Out] uint8_t% year,`
`[System::Runtime::InteropServices::Out] uint8_t% month,`
`[System::Runtime::InteropServices::Out] uint8_t% day,`
`[System::Runtime::InteropServices::Out] uint8_t% hour,`
`[System::Runtime::InteropServices::Out] uint8_t% minute,`
`[System::Runtime::InteropServices::Out] uint8_t% second)`

Get the RTC

Parameters

<i>year</i>	the year
<i>month</i>	the month
<i>day</i>	the day
<i>hour</i>	the hour
<i>minute</i>	the minute
<i>second</i>	the second

11.85.3.10 GetStateDebugData() `uint32_t GetStateDebugData (`
`uint16_t coil)`

get the debug queue state

Parameters

<i>coil</i>	the coil
-------------	----------

Returns

the state

11.85.3.11 GetStateEventData() `uint32_t GetStateEventData ()`

get the event queue state

Returns

the state

11.85.3.12 SetImplantCurrentSetpoint() `void SetImplantCurrentSetpoint (`
`uint16_t coil,`
`uint32_t current)`

gets the implant current setpoint

Parameters

<i>coil</i>	the coil
<i>current</i>	the current

11.85.3.13 SetPowerStrength() `void SetPowerStrength (`
`uint16_t coil,`
`uint32_t power)`

gets the power for the trigger pulses

Parameters

<i>coil</i>	the coil
<i>power</i>	the power

11.85.3.14 SetRTC() `void SetRTC (`
`uint8_t year,`
`uint8_t month,`
`uint8_t day,`
`uint8_t hour,`
`uint8_t minute,`
`uint8_t second)`

Set the RTC

Parameters

<i>year</i>	the year
<i>month</i>	the month
<i>day</i>	the day
<i>hour</i>	the hour
<i>minute</i>	the minute
<i>second</i>	the second

11.85.3.15 SetStateDebugData() void SetStateDebugData (
 uint16_t *coil*,
 uint32_t *state*)

clears/starts/stops the debug queue for a certain coil

Parameters

<i>coil</i>	the coil
<i>state</i>	clear/start/stop

11.85.3.16 SetStateEventData() void SetStateEventData (
 uint32_t *state*)

clears/starts/stops the event queue for a certain coil

Parameters

<i>state</i>	clear/start/stop
--------------	------------------

11.85.3.17 SwitchOnOff() void SwitchOnOff (
 uint16_t *coil*,
 uint32_t *on*)

switched the coild on of

Parameters

<i>coil</i>	the coil
<i>on</i>	0 = off, 1 = on


11.85.4 Property Documentation

11.85.4.1 RFFunction `CRFFFunctionNet^ RFFunction [get]`

11.86 CPositionImpDeviceNet Class Reference

`CPositionImpDeviceNet` is the class to access the Position/Imp devices

Inheritance diagram for `CPositionImpDeviceNet`:

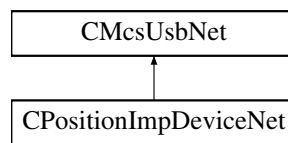

Public Member Functions

- `CPositionImpDeviceNet ()`
Initializes a new instance of the `CPositionImpDeviceNet` class.
- virtual `~CPositionImpDeviceNet ()`
- `!CPositionImpDeviceNet ()`
- void `ConnectImp (uint32_t id)`
Connect to a Imp device with a certain ID
- uint32_t `ConnectedImp ()`
The ID of the connected Imp device
- int32_t `GetRFFrequency ()`
Gets currently used RF frequency
- void `SetRFFrequency (int32_t frequency)`
Sets the current RF frequency
- uint32_t `GetDeviceList (int32_t index)`
Gets the device list
- void `SetDeviceList (int32_t index, uint32_t id)`
Sets the device list
- uint32_t `GetImpId ()`
Gets the ID of the impedance measure device
- void `SetImpId (uint32_t id)`
Sets the ID of the impedance measure device

Additional Inherited Members

11.86.1 Detailed Description

`CPositionImpDeviceNet` is the class to access the Position/Imp devices

11.86.2 Constructor & Destructor Documentation

11.86.2.1 CPositionImpDeviceNet() `CPositionImpDeviceNet ()`

Initializes a new instance of the `CPositionImpDeviceNet` class.

11.86.2.2 ~CPositionImpDeviceNet() `virtual ~CPositionImpDeviceNet () [virtual]`

11.86.2.3 "!CPositionImpDeviceNet() `!CPositionImpDeviceNet ()`

11.86.3 Member Function Documentation

11.86.3.1 ConnectedImp() `uint32_t ConnectedImp ()`

The ID of the connected Imp device

Returns

The ID

11.86.3.2 ConnectImp() `void ConnectImp (uint32_t id)`

Connect to a Imp device with a certain ID

Parameters

<i>id</i>	The ID
-----------	--------

11.86.3.3 GetDeviceList() `uint32_t GetDeviceList (int32_t index)`

Gets the device list

Parameters

<i>index</i>	the index
--------------	-----------

Returns

the ID

11.86.3.4 GetImpId() `uint32_t GetImpId ()`

Gets the ID of the impedance measure device

Returns

the ID

11.86.3.5 GetRFFrequency() `int32_t GetRFFrequency ()`

Gets currently used RF frequency

Returns

The frequency

11.86.3.6 SetDeviceList() `void SetDeviceList (`
 `int32_t index,`
 `uint32_t id)`

Sets the device list

Parameters

<i>index</i>	the index
<i>id</i>	the ID

11.86.3.7 SetImpId() `void SetImpId (`
 `uint32_t id)`

Sets the ID of the impedance measure device

Parameters

<i>id</i>	the ID
-----------	--------

11.86.3.8 SetRFFrequency() `void SetRFFrequency (`
`int32_t frequency)`


Sets the current RF frequency

Parameters

<i>frequency</i>	The frequency
------------------	---------------

11.87 CPPCDeviceNet Class Reference

Inheritance diagram for CPPCDeviceNet:

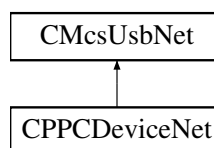

Public Member Functions

- [CPPCDeviceNet](#) (void)

Properties

- [CPPCFunctionNet](#)[^] [PPCFunction](#) [get]
- [CMcsBusNet](#)[^] [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]
- [CMcsBus_SensorNet](#)[^] [McsBus_Sensor](#) [get]

Additional Inherited Members

11.87.1 Constructor & Destructor Documentation

11.87.1.1 CPPCDeviceNet() `CPPCDeviceNet (`
`void)`

11.87.2 Property Documentation

11.87.2.1 McsBus [CMcsBusNet](#)[^] McsBus [get]

11.87.2.2 McsBus_MotorControl [CMcsBus_MotorControlNet](#)[^] McsBus_MotorControl [get]


11.87.2.3 McsBus_Sensor [CMcsBus_SensorNet](#)[^] McsBus_Sensor [get]

11.87.2.4 PPCFunction [CPPPCFunctionNet](#)[^] PPCFunction [get]

11.88 CPPCFunctionNet Class Reference

[CPPCFunctionNet](#) is the class to access the PPC (high precision Patch Peristaltic patch Pump

Inheritance diagram for CPPCFunctionNet:

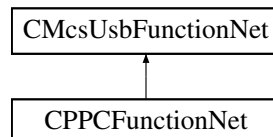

Public Member Functions

- [CPPCFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pPPCFunctionPointer↔ Container)
Initializes a new instance of the [CPPCFunctionNet](#) class.
- [CPPCFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CPPCFunctionNet](#) ()
- [!CPPCFunctionNet](#) ()
- int [GetPumpSpeedUnit](#) (uint16_t channel)
Reads the Pump Speed Unit
- void [SetPumpSpeedUnit](#) (uint16_t channel, int SpeedUnit)
Writes the Pump Speed Unit
- [PP_Pump_Mode_Type_EnumNet](#) [GetPumpModeType](#) (uint16_t channel)
Reads the Pump Mode Type.
- void [SetPumpModeType](#) (uint16_t channel, [PP_Pump_Mode_Type_EnumNet](#) PumpMode)
Writes the config string from the device.

- void [GetAnalogVoltageRange](#) (uint16_t channel, [System::Runtime::InteropServices::Out]uint16_t% min_voltage, [System::Runtime::InteropServices::Out]uint16_t% max_voltage)
Reads the Analog Input Voltage Range
- void [SetAnalogVoltageRange](#) (uint16_t channel, uint16_t min_voltage, uint16_t max_voltage)
Writes the Analog Input Voltage Range
- void [GetPressureRange](#) (uint16_t channel, [System::Runtime::InteropServices::Out]int32_t% lower_pressure, [System::Runtime::InteropServices::Out]int32_t% upper_pressure)
Get the pressure range that is used between the analog voltage or the digital states
- void [SetPressureRange](#) (uint16_t channel, int32_t lower_pressure, int32_t upper_pressure)
Get the pressure range that is used between the analog voltage or the digital states
- uint16_t [GetSupplyVoltage](#) ()
Reads the current supply voltage in mV
- uint16_t [GetAnalogVoltage](#) (uint16_t channel)
Reads the current analog voltage
- uint16_t [GetDigitalIn](#) (uint16_t channel)
Reads the digital input state
- int [GetValveActive](#) (uint16_t valve)
Gets the valve active/inactive state
- void [SetValveActive](#) (uint16_t valve, int valveActive)
Sets the valve active/inactive state
- void [SetPressureOffset](#) ()
Sets the pressure offset
- void [LoadPressure](#) (int32_t pressure, uint32_t options)
Loads the reservoir with a pressure
- void [IsBusy](#) ([System::Runtime::InteropServices::Out]int16_t% task, [System::Runtime::InteropServices::Out]int16_t% wait)
Is the PPC busy with a task
- void [FirePressurePulse](#) (int32_t duration, int32_t nextpressure)
Fire a pressure pulse from the reservoir
- int32_t [MeasureReservoir](#) ()
Measures the reservoir pressure

Additional Inherited Members

11.88.1 Detailed Description

[CPPCFunctionNet](#) is the class to access the PPC (high precision Patch Peristaltic patch Pump

11.88.2 Constructor & Destructor Documentation

11.88.2.1 CPPCFunctionNet() [1/2] [CPPCFunctionNet](#) (
 CMcsUsbNet^ mcsusb,
 CMcsUsbFunctionPointerContainer^ pPPCFunctionPointerContainer)

Initializes a new instance of the [CPPCFunctionNet](#) class.

11.88.2.2 CPPCFunctionNet() [2/2] `CPPCFunctionNet (`
`CMcsUsbNet^ mcsusb)`

11.88.2.3 ~CPPCFunctionNet() `virtual ~CPPCFunctionNet () [virtual]`

11.88.2.4 "!CPPCFunctionNet() `!CPPCFunctionNet ()`

11.88.3 Member Function Documentation

11.88.3.1 FirePressurePulse() `void FirePressurePulse (`
`int32_t duration,`
`int32_t nextpressure)`

Fire a pressure pulse from the reservoir

Parameters

<i>duration</i>	The pulse duration (valves open)
<i>nextpressure</i>	The next pressure

11.88.3.2 GetAnalogVoltage() `uint16_t GetAnalogVoltage (`
`uint16_t channel)`

Reads the current analog voltage

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Analog Voltage

11.88.3.3 GetAnalogVoltageRange() `void GetAnalogVoltageRange (`
`uint16_t channel,`

```
[System::Runtime::InteropServices::Out] uint16_t% min_voltage,  
[System::Runtime::InteropServices::Out] uint16_t% max_voltage )
```

Reads the Analog Input Voltage Range

Parameters

<i>channel</i>	The Channel Number
<i>min_voltage</i>	The voltage that should be seen as the minimum voltage
<i>max_voltage</i>	The voltage that should be seen as the maximum voltage

11.88.3.4 GetDigitalIn() `uint16_t GetDigitalIn (`
`uint16_t channel)`

Reads the digital input state

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Digital State

11.88.3.5 GetPressureRange() `void GetPressureRange (`
`uint16_t channel,`
`[System::Runtime::InteropServices::Out] int32_t% lower_pressure,`
`[System::Runtime::InteropServices::Out] int32_t% upper_pressure)`

Get the pressure range that is used between the analog voltage or the digital states

Parameters

<i>channel</i>	The Channel Number
<i>lower_pressure</i>	The lower border of the pressure range
<i>upper_pressure</i>	The upper border of the pressure range

11.88.3.6 GetPumpModeType() `PP_Pump_Mode_Type_EnumNet GetPumpModeType (`
`uint16_t channel)`

Reads the Pump Mode Type.

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Pump Mode Type.

11.88.3.7 GetPumpSpeedUnit() `int GetPumpSpeedUnit (uint16_t channel)`

Reads the Pump Speed Unit

Parameters

<i>channel</i>	The Channel Number
----------------	--------------------

Returns

The Speed Unit

11.88.3.8 GetSupplyVoltage() `uint16_t GetSupplyVoltage ()`

Reads the current supply voltage in mV

Returns

The supply voltage

11.88.3.9 GetValveActive() `int GetValveActive (uint16_t valve)`

Gets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve state

11.88.3.10 IsBusy() `void IsBusy ([System::Runtime::InteropServices::Out] int16_t% task, [System::Runtime::InteropServices::Out] int16_t% wait)`

Is the PPC busy with a task

Parameters

<i>task</i>	The task state
<i>wait</i>	The wait state

11.88.3.11 LoadPressure() `void LoadPressure (`
 `int32_t pressure,`
 `uint32_t options)`

Loads the reservoir with a pressure

Parameters

<i>pressure</i>	The pressure
<i>options</i>	The options: end with 0=regulate on patch 1=regulate on reservoir

11.88.3.12 MeasureReservoir() `int32_t MeasureReservoir ()`

Measures the reservoir pressure

Returns

The pressure

11.88.3.13 SetAnalogVoltageRange() `void SetAnalogVoltageRange (`
 `uint16_t channel,`
 `uint16_t min_voltage,`
 `uint16_t max_voltage)`

Writes the Analog Input Voltage Range

Parameters

<i>channel</i>	The Channel Number
<i>min_voltage</i>	The voltage that should be seen as the minimum voltage
<i>max_voltage</i>	The voltage that should be seen as the maximum voltage

11.88.3.14 SetPressureOffset() `void SetPressureOffset ()`

Sets the pressure offset

11.88.3.15 SetPressureRange() `void SetPressureRange (`
 `uint16_t channel,`
 `int32_t lower_pressure,`
 `int32_t upper_pressure)`

Get the pressure range that is used between the analog voltage or the digital states

Parameters

<i>channel</i>	The Channel Number
<i>lower_pressure</i>	The lower border of the pressure range
<i>upper_pressure</i>	The upper border of the pressure range

11.88.3.16 SetPumpModeType() `void SetPumpModeType (`
 `uint16_t channel,`
 `PP_Pump_Mode_Type_EnumNet PumpMode)`

Writes the config string from the device.

Parameters

<i>channel</i>	The Channel Number
<i>PumpMode</i>	The Pump Mode Type.

11.88.3.17 SetPumpSpeedUnit() `void SetPumpSpeedUnit (`
 `uint16_t channel,`
 `int SpeedUnit)`

Writes the Pump Speed Unit

Parameters

<i>channel</i>	The Channel Number
<i>SpeedUnit</i>	The Speed Unit

11.88.3.18 SetValveActive() `void SetValveActive (`
 `uint16_t valve,`
 `int valveActive)`


Sets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
<i>valveActive</i>	The valve state

11.89 CPPS_DeviceNet Class Reference

Inheritance diagram for CPPS_DeviceNet:

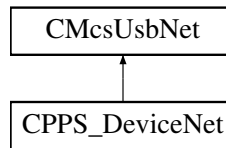

Public Member Functions

- [CPPS_DeviceNet](#) (void)

Properties

- [CPPS_FunctionNet](#)[^] [PPS_Function](#) [get]
- [CMcsBusNet](#)[^] [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]
- [CMcsBus_SensorNet](#)[^] [McsBus_Sensor](#) [get]

Additional Inherited Members

11.89.1 Constructor & Destructor Documentation

11.89.1.1 [CPPS_DeviceNet\(\)](#) [CPPS_DeviceNet](#) (
void)

11.89.2 Property Documentation

11.89.2.1 **McsBus** [CMcsBusNet](#)[^] [McsBus](#) [get]


11.89.2.2 **McsBus_MotorControl** [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]

11.89.2.3 **McsBus_Sensor** [CMcsBus_SensorNet](#)[^] [McsBus_Sensor](#) [get]

11.89.2.4 PPS_Function CPPS_FunctionNet^ PPS_Function [get]

11.90 CPPS_FunctionNet Class Reference

Inheritance diagram for CPPS_FunctionNet:

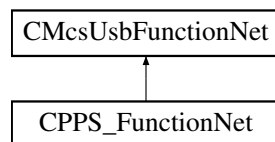

Public Member Functions

- [CPPS_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ cPPS_FunctionPointerContainer)
- [CPPS_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- void [SetPumpMaxSpeed](#) (unsigned short index, unsigned short maxspeed)
- unsigned short [GetPumpMaxSpeed](#) (unsigned short index)
- void [SetPumpSpeedUnit](#) (unsigned short index, int speedunit)
- int [GetPumpSpeedUnit](#) (unsigned short index)
- void [SetPumpModeType](#) (unsigned short index, [PP_Pump_Mode_Type_EnumNet](#) type)
- [PP_Pump_Mode_Type_EnumNet](#) [GetPumpModeType](#) (unsigned short index)
- void [SetPumpCouple](#) (unsigned int i)
- unsigned int [GetPumpCouple](#) ()
- void [SetPumpEnableSpeedRatio](#) (unsigned int enable)
- unsigned int [GetPumpEnableSpeedRatio](#) ()
- void [SetPumpManualOnOff](#) (unsigned short index, unsigned int onoff)
- unsigned int [GetPumpManualOnOff](#) (unsigned short index)
- void [SetPumpFunctionSpeeds](#) (unsigned short index, short offspeed, short onspeed)
- void [GetPumpFunctionSpeeds](#) (unsigned short index, [System::Runtime::InteropServices::Out]short% offspeed, [System::Runtime::InteropServices::Out]short% onspeed)
- void [SetPumpSpeedRatio](#) (int ratio)
- int [GetPumpSpeedRatio](#) ()
- void [SetPumpFastOnOff](#) (unsigned short index, unsigned int onoff)
- unsigned int [GetPumpFastOnOff](#) (unsigned short index)
- void [SetPumpFastSpeed](#) (unsigned short index, short fastspeed)
- short [GetPumpFastSpeed](#) (unsigned short index)
- void [SetAnalogVoltages](#) (unsigned short index, unsigned short minvoltage, unsigned short maxvoltage)
- void [GetAnalogVoltages](#) (unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% minvoltage, [System::Runtime::InteropServices::Out]unsigned short% maxvoltage)
- void [SetUseBubble](#) (unsigned short index, unsigned int usebubble)
- unsigned int [GetUseBubble](#) (unsigned short index)
- unsigned short [GetSupplyVoltage](#) ()
- unsigned short [GetAnalogVoltage](#) (unsigned short index)
- unsigned short [GetDigitalIn](#) (unsigned short index)
- unsigned short [GetBubbleState](#) ()

Additional Inherited Members

11.90.1 Constructor & Destructor Documentation

11.90.1.1 CPPS_FunctionNet() [1/2] `CPPS_FunctionNet (`
 `CMcsUsbNet^ mcsusb,`
 `CMcsUsbFunctionPointerContainer^ cPPS_FunctionPointerContainer)`

11.90.1.2 CPPS_FunctionNet() [2/2] `CPPS_FunctionNet (`
 `CMcsUsbNet^ mcsusb)`

11.90.2 Member Function Documentation

11.90.2.1 GetAnalogVoltage() `unsigned short GetAnalogVoltage (`
 `unsigned short index)`

11.90.2.2 GetAnalogVoltages() `void GetAnalogVoltages (`
 `unsigned short index,`
 `[System::Runtime::InteropServices::Out] unsigned short% minvoltage,`
 `[System::Runtime::InteropServices::Out] unsigned short% maxvoltage)`

11.90.2.3 GetBubbleState() `unsigned short GetBubbleState ()`

11.90.2.4 GetDigitalIn() `unsigned short GetDigitalIn (`
 `unsigned short index)`

11.90.2.5 GetPumpCouple() `unsigned int GetPumpCouple ()`

11.90.2.6 GetPumpEnableSpeedRatio() `unsigned int GetPumpEnableSpeedRatio ()`

11.90.2.7 GetPumpFastOnOff() `unsigned int GetPumpFastOnOff (`
 `unsigned short index)`

11.90.2.8 GetPumpFastSpeed() short GetPumpFastSpeed (
 unsigned short *index*)

11.90.2.9 GetPumpFunctionSpeeds() void GetPumpFunctionSpeeds (
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] short% *offspeed*,
 [System::Runtime::InteropServices::Out] short% *onspeed*)

11.90.2.10 GetPumpManualOnOff() unsigned int GetPumpManualOnOff (
 unsigned short *index*)

11.90.2.11 GetPumpMaxSpeed() unsigned short GetPumpMaxSpeed (
 unsigned short *index*)

11.90.2.12 GetPumpModeType() [PP_Pump_Mode_Type_EnumNet](#) GetPumpModeType (
 unsigned short *index*)

11.90.2.13 GetPumpSpeedRatio() int GetPumpSpeedRatio ()

11.90.2.14 GetPumpSpeedUnit() int GetPumpSpeedUnit (
 unsigned short *index*)

11.90.2.15 GetSupplyVoltage() unsigned short GetSupplyVoltage ()

11.90.2.16 GetUseBubble() unsigned int GetUseBubble (
 unsigned short *index*)

11.90.2.17 SetAnalogVoltages() void SetAnalogVoltages (
 unsigned short *index*,
 unsigned short *minvoltage*,
 unsigned short *maxvoltage*)

11.90.2.18 SetPumpCouple() void SetPumpCouple (
 unsigned int *i*)

11.90.2.19 SetPumpEnableSpeedRatio() void SetPumpEnableSpeedRatio (
 unsigned int *enable*)

11.90.2.20 SetPumpFastOnOff() void SetPumpFastOnOff (
 unsigned short *index*,
 unsigned int *onoff*)

11.90.2.21 SetPumpFastSpeed() void SetPumpFastSpeed (
 unsigned short *index*,
 short *fastspeed*)

11.90.2.22 SetPumpFunctionSpeeds() void SetPumpFunctionSpeeds (
 unsigned short *index*,
 short *offspeed*,
 short *onspeed*)

11.90.2.23 SetPumpManualOnOff() void SetPumpManualOnOff (
 unsigned short *index*,
 unsigned int *onoff*)

11.90.2.24 SetPumpMaxSpeed() void SetPumpMaxSpeed (
 unsigned short *index*,
 unsigned short *maxspeed*)

11.90.2.25 SetPumpModeType() `void SetPumpModeType (`
 `unsigned short index,`
 `PP_Pump_Mode_Type_EnumNet type)`

11.90.2.26 SetPumpSpeedRatio() `void SetPumpSpeedRatio (`
 `int ratio)`


11.90.2.27 SetPumpSpeedUnit() `void SetPumpSpeedUnit (`
 `unsigned short index,`
 `int speedunit)`

11.90.2.28 SetUseBubble() `void SetUseBubble (`
 `unsigned short index,`
 `unsigned int usebubble)`

11.91 CPPSDeviceNet Class Reference

CPPS4plus1DeviceNet is the to control the MCS HLA device

Inheritance diagram for CPPSDeviceNet:

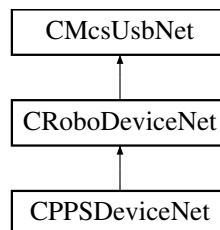

Public Member Functions

- [CPPSDeviceNet](#) (void)

Additional Inherited Members

11.91.1 Detailed Description

CPPS4plus1DeviceNet is the to control the MCS HLA device


11.91.2 Constructor & Destructor Documentation

11.91.2.1 CPPSDeviceNet() `CPPSDeviceNet (void)`

11.92 CProgramPressureCurveNet Class Reference

`CProgramPressureCurveNet` is the class to program pressure curves

Inheritance diagram for `CProgramPressureCurveNet`:

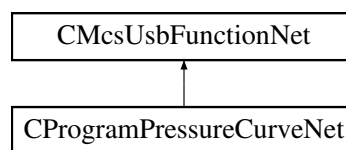

Public Member Functions

- `CProgramPressureCurveNet (CMcsUsbNet^ mcsusb)`
Initializes a new instance of the `CPulseGeneratorFunctionNet` class.
- virtual `~CProgramPressureCurveNet (void)`
- `!CProgramPressureCurveNet (void)`
- void `Program` (unsigned char busnumber, unsigned char busaddress, int32_t channel, array< int32_t >^ pressures, array< int32_t >^ steps, array< int16_t >^ durations)
- void `SetRepeats` (unsigned char busnumber, unsigned char busaddress, int32_t channel, uint32_t repeats)
- unsigned int `GetRepeats` (unsigned char busnumber, unsigned char busaddress, int32_t channel)

Additional Inherited Members

11.92.1 Detailed Description

`CProgramPressureCurveNet` is the class to program pressure curves

11.92.2 Constructor & Destructor Documentation

11.92.2.1 CProgramPressureCurveNet() `CProgramPressureCurveNet (CMcsUsbNet^ mcsusb)`

Initializes a new instance of the `CPulseGeneratorFunctionNet` class.

11.92.2.2 `~CProgramPressureCurveNet()` `virtual ~CProgramPressureCurveNet (`
`void) [virtual]`

11.92.2.3 `!"CProgramPressureCurveNet()` `!CProgramPressureCurveNet (`
`void)`

11.92.3 Member Function Documentation

11.92.3.1 `GetRepeats()` `unsigned int GetRepeats (`
`unsigned char busnumber,`
`unsigned char busaddress,`
`int32_t channel)`


11.92.3.2 `Program()` `void Program (`
`unsigned char busnumber,`
`unsigned char busaddress,`
`int32_t channel,`
`array< int32_t >^ pressures,`
`array< int32_t >^ steps,`
`array< int16_t >^ durations)`

11.92.3.3 `SetRepeats()` `void SetRepeats (`
`unsigned char busnumber,`
`unsigned char busaddress,`
`int32_t channel,`
`uint32_t repeats)`

11.93 CPulseGeneratorFunctionNet Class Reference

[CPulseGeneratorFunctionNet](#) is the class to control the pulse generator for video tracking

Inheritance diagram for [CPulseGeneratorFunctionNet](#):

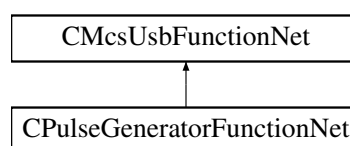

Public Member Functions

- [CPulseGeneratorFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pPulseGeneratorFunctionPointerContainer)
Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.
- [CPulseGeneratorFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual [~CPulseGeneratorFunctionNet](#) ()
- [ICPulseGeneratorFunctionNet](#) ()
- [int32_t](#) [GetPeriod](#) ([int32_t](#) generator_number)
Reads the generator period
- void [SetPeriod](#) ([int32_t](#) generator_number, [int32_t](#) period_in_samples)
Writes the generator period
- [int32_t](#) [GetPulseLength](#) ([int32_t](#) generator_number)
Reads the generator pulse length
- void [SetPulseLength](#) ([int32_t](#) generator_number, [int32_t](#) pulselength_in_10us)
Writes the generator pulse length
- void [GetModeSelect](#) ([int32_t](#) generator_number, [System::Runtime::InteropServices::Out][PulseGenerator_Mode_EnumNet](#)% mode, [System::Runtime::InteropServices::Out][int32_t](#)% digitalchannel)
Reads the generator mode
- void [SetModeSelect](#) ([int32_t](#) generator_number, [PulseGenerator_Mode_EnumNet](#) mode, [int32_t](#) digitalchannel)
Writes the generator mode

Additional Inherited Members

11.93.1 Detailed Description

[CPulseGeneratorFunctionNet](#) is the class to control the pulse generator for video tracking

11.93.2 Constructor & Destructor Documentation

- 11.93.2.1 [CPulseGeneratorFunctionNet](#)() [1/2]** [CPulseGeneratorFunctionNet](#) (
[CMcsUsbNet](#)^ mcsusb,
[CMcsUsbFunctionPointerContainer](#)^ pPulseGeneratorFunctionPointerContainer)

Initializes a new instance of the [CPulseGeneratorFunctionNet](#) class.

- 11.93.2.2 [CPulseGeneratorFunctionNet](#)() [2/2]** [CPulseGeneratorFunctionNet](#) (
[CMcsUsbNet](#)^ mcsusb)

- 11.93.2.3 [~CPulseGeneratorFunctionNet](#)()** virtual [~CPulseGeneratorFunctionNet](#) () [virtual]

11.93.2.4 `!CPulseGeneratorFunctionNet()` `!CPulseGeneratorFunctionNet ()`

11.93.3 Member Function Documentation

11.93.3.1 GetModeSelect() `void GetModeSelect (`
 `int32_t generator_number,`
 `[System::Runtime::InteropServices::Out] PulseGenerator_Mode_EnumNet% mode,`
 `[System::Runtime::InteropServices::Out] int32_t% digitalchannel)`

Reads the generator mode

Parameters

<i>generator_number</i>	The generator number
<i>mode</i>	The generator mode
<i>digitalchannel</i>	The digital in channel used as gate

11.93.3.2 GetPeriod() `int32_t GetPeriod (`
 `int32_t generator_number)`

Reads the generator period

Parameters

<i>generator_number</i>	The generator number
-------------------------	----------------------

Returns

The period

11.93.3.3 GetPulseLength() `int32_t GetPulseLength (`
 `int32_t generator_number)`

Reads the generator pulse length

Parameters

<i>generator_number</i>	The generator number
-------------------------	----------------------

Returns

The pulse length

11.93.3.4 SetModeSelect() `void SetModeSelect (`
`int32_t generator_number,`
`PulseGenerator_Mode_EnumNet mode,`
`int32_t digitalchannel)`

Writes the generator mode

Parameters

<i>generator_number</i>	The generator number
<i>mode</i>	The generator mode
<i>digitalchannel</i>	The digital in channel used as gate

11.93.3.5 SetPeriod() `void SetPeriod (`
`int32_t generator_number,`
`int32_t period_in_samples)`

Writes the generator period

Parameters

<i>generator_number</i>	The generator number
<i>period_in_samples</i>	The period

11.93.3.6 SetPulseLength() `void SetPulseLength (`
`int32_t generator_number,`
`int32_t pulselength_in_10us)`


Writes the generator pulse length

Parameters

<i>generator_number</i>	The generator number
<i>pulselength_in_10us</i>	The pulse length

11.94 CRadioControlledDevicesNet Class Reference

Inheritance diagram for CRadioControlledDevicesNet:

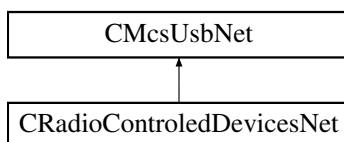

Public Member Functions

- [CRadioControlledDevicesNet](#) (void)
- bool [HasRadioControl](#) ()
- array< unsigned short > ^ [GetDeviceNames](#) ()
- void [ConnectDevice](#) (unsigned short sn)
- void [DisConnectDevice](#) ()
- bool [StillConnected](#) ()
- void [SetFrequency](#) (unsigned short frequency)
- unsigned short [GetFrequency](#) ()

Protected Member Functions

- [CRadioControlledDevicesNet](#) (CRadioControlledDevices *pRadioControlled)

Additional Inherited Members

11.94.1 Constructor & Destructor Documentation

11.94.1.1 CRadioControlledDevicesNet() [1/2] [CRadioControlledDevicesNet](#) (
void)

11.94.1.2 CRadioControlledDevicesNet() [2/2] [CRadioControlledDevicesNet](#) (
CRadioControlledDevices * *pRadioControlled*) [protected]

11.94.2 Member Function Documentation

11.94.2.1 ConnectDevice() void [ConnectDevice](#) (
unsigned short *sn*)

11.94.2.2 DisConnectDevice() void [DisConnectDevice](#) ()

11.94.2.3 GetDeviceNames() `array<unsigned short> ^ GetDeviceNames ()`

11.94.2.4 GetFrequency() `unsigned short GetFrequency ()`

11.94.2.5 HasRadioControl() `bool HasRadioControl ()`

11.94.2.6 SetFrequency() `void SetFrequency (`
`unsigned short frequency)`

11.94.2.7 StillConnected() `bool StillConnected ()`

11.95 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference

Public Member Functions

- [CRegionOfInterestRect](#) (int left, int top, int right, int bottom)
- [CRegionOfInterestRect](#) ^ [DeepCopy](#) ()

Public Attributes

- int [m_Left](#)
- int [m_Top](#)
- int [m_Right](#)
- int [m_Bottom](#)

11.95.1 Constructor & Destructor Documentation

11.95.1.1 CRegionOfInterestRect() `CRegionOfInterestRect (`
`int left,`
`int top,`
`int right,`
`int bottom)`

11.95.2 Member Function Documentation

11.95.2.1 DeepCopy() [CRegionOfInterestRect](#) ^ DeepCopy ()

11.95.3 Member Data Documentation

11.95.3.1 m_Bottom `int m_Bottom`


11.95.3.2 m_Left `int m_Left`

11.95.3.3 m_Right `int m_Right`

11.95.3.4 m_Top `int m_Top`

11.96 CRetinaLedDeviceNet Class Reference

Inheritance diagram for CRetinaLedDeviceNet:

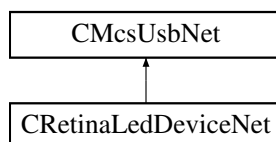

Public Member Functions

- [CRetinaLedDeviceNet](#) ()
- [~CRetinaLedDeviceNet](#) ()
- unsigned int [SetTrigger](#) (int enable)
- unsigned int [SetLED](#) (unsigned long long pattern)
- unsigned int [SetTablepointer](#) (int position)
- unsigned int [GetTablepointer](#) (int % position)
- unsigned int [ClearTable](#) ()
- unsigned int [AddTableEntry](#) (unsigned long long pattern)
- unsigned int [AddLoopEntry](#) (unsigned short repeats, unsigned short steps_back)
- unsigned int [SetRepeat](#) (int repeat)
- unsigned int [SetLumi](#) (int lumi)
- unsigned int [SetPersistency](#) (unsigned int persistency)

Additional Inherited Members

11.96.1 Constructor & Destructor Documentation

11.96.1.1 CRetinaLedDeviceNet() `CRetinaLedDeviceNet ()`

11.96.1.2 ~CRetinaLedDeviceNet() `~CRetinaLedDeviceNet ()`

11.96.2 Member Function Documentation

11.96.2.1 AddLoopEntry() `unsigned int AddLoopEntry (unsigned short repeats, unsigned short steps_back)`

11.96.2.2 AddTableEntry() `unsigned int AddTableEntry (unsigned long long pattern)`

11.96.2.3 ClearTable() `unsigned int ClearTable ()`

11.96.2.4 GetTablepointer() `unsigned int GetTablepointer (int % position)`

11.96.2.5 SetLED() `unsigned int SetLED (unsigned long long pattern)`

11.96.2.6 SetLumi() `unsigned int SetLumi (int lumi)`

11.96.2.7 SetPersistence() unsigned int SetPersistence (
 unsigned int *persistence*)

11.96.2.8 SetRepeat() unsigned int SetRepeat (
 int *repeat*)


11.96.2.9 SetTablepointer() unsigned int SetTablepointer (
 int *position*)

11.96.2.10 SetTrigger() unsigned int SetTrigger (
 int *enable*)

11.97 CRFFunctionNet Class Reference

[CRFFunctionNet](#) is the class to control RF devices

Inheritance diagram for CRFFunctionNet:

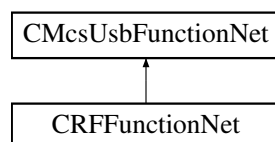

Public Member Functions

- [CRFFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pRFFunctionPointerContainer)
 Initializes a new instance of the [CRFFunctionNet](#) class.
- [CRFFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CRFFunctionNet](#) ()
- [!CRFFunctionNet](#) ()
- uint32_t [GetBaseFrequency](#) ([CFirmwareDestinationNet](#) destination)
 gets the base advertise frequency
- void [SetBaseFrequency](#) ([CFirmwareDestinationNet](#) destination, uint32_t frequency)
 sets the base advertise frequency
- uint32_t [GetWorkingFrequency](#) ()
 gets the working frequency
- void [SetWorkingFrequency](#) (uint32_t frequency)
 sets the working frequency
- array< uint32_t > [^] [GetAvailableDeviceListEx](#) (int list_Length)
 get a list of available devices
- array< uint32_t > [^] [GetAvailableDeviceList](#) ()

- get a list of available devices*
- array< uint32_t > ^ [GetAvailableStateListEx](#) (int list_Length)
- get a list of the states of the available devices*
- array< uint32_t > ^ [GetAvailableStateList](#) ()
- get a list of the states of the available devices*
- void [Connect](#) (uint32_t sn)
- connect to a RF device, use 0 to disconnect*
- uint32_t [GetConnectedDevice](#) ()
- get connect RF device, 0 = no device connected*
- uint32_t [GetState](#) ()
- get connection state*
- void [SetTestMode](#) (uint32_t mode)
- set test mode*
- uint32_t [GetTestMode](#) ()
- gets test mode*

Additional Inherited Members

11.97.1 Detailed Description

[CRFFunctionNet](#) is the class to control RF devices

11.97.2 Constructor & Destructor Documentation

11.97.2.1 [CRFFunctionNet\(\)](#) [1/2] [CRFFunctionNet](#) (
 CMcsUsbNet^ mcsusb,
 CMcsUsbFunctionPointerContainer^ pRFFunctionPointerContainer)

Initializes a new instance of the [CRFFunctionNet](#) class.

11.97.2.2 [CRFFunctionNet\(\)](#) [2/2] [CRFFunctionNet](#) (
 CMcsUsbNet^ mcsusb)

11.97.2.3 [~CRFFunctionNet\(\)](#) virtual [~CRFFunctionNet](#) () [virtual]

11.97.2.4 ["!CRFFunctionNet\(\)](#) [!CRFFunctionNet](#) ()

11.97.3 Member Function Documentation

11.97.3.1 [Connect\(\)](#) void [Connect](#) (
 uint32_t sn)

connect to a RF device, use 0 to disconnect

Parameters

<i>sn</i>	the serial number
-----------	-------------------

11.97.3.2 GetAvailableDeviceList() `array<uint32_t> ^ GetAvailableDeviceList ()`

get a list of available devices

Returns

array of devices

11.97.3.3 GetAvailableDeviceListEx() `array<uint32_t> ^ GetAvailableDeviceListEx (
int list_Length)`

get a list of available devices

Parameters

<i>list_Length</i>	The maximal length of list.
--------------------	-----------------------------

Returns

array of devices

11.97.3.4 GetAvailableStateList() `array<uint32_t> ^ GetAvailableStateList ()`

get a list of the states of the available devices

Returns

array of states

11.97.3.5 GetAvailableStateListEx() `array<uint32_t> ^ GetAvailableStateListEx (
int list_Length)`

get a list of the states of the available devices

Parameters

<i>list_Length</i>	The maximal length of list.
--------------------	-----------------------------

Returns

array of states

11.97.3.6 GetBaseFrequency() `uint32_t GetBaseFrequency (CFirmwareDestinationNet destination)`

gets the base advertise frequency

Parameters

<i>destination</i>	the destination to query
--------------------	--------------------------

Returns

the frequency

11.97.3.7 GetConnectedDevice() `uint32_t GetConnectedDevice ()`

get connect RF device, 0 = no device connected

Returns

the serial number

11.97.3.8 GetState() `uint32_t GetState ()`

get connection state

Returns

the state

11.97.3.9 GetTestMode() `uint32_t GetTestMode ()`

gets test mode

Returns

the mode

11.97.3.10 GetWorkingFrequency() `uint32_t GetWorkingFrequency ()`

gets the working frequency

Returns

the frequency

11.97.3.11 SetBaseFrequency() `void SetBaseFrequency (
 CFirmwareDestinationNet destination,
 uint32_t frequency)`

sets the base advertise frequency

Parameters

<i>destination</i>	the destination to set
<i>frequency</i>	the frequency

11.97.3.12 SetTestMode() `void SetTestMode (
 uint32_t mode)`

set test mode

Parameters

<i>mode</i>	the mode
-------------	----------

11.97.3.13 SetWorkingFrequency() `void SetWorkingFrequency (
 uint32_t frequency)`


sets the working frequency

Parameters

<i>frequency</i>	the frequency
------------------	---------------

11.98 CRobo_FYIPProgram_FunctionNet Class Reference

Inheritance diagram for CRobo_FYIPProgram_FunctionNet:

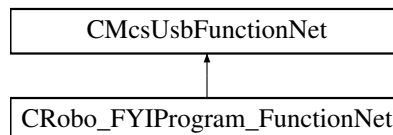

Public Member Functions

- [CRobo_FYIPProgram_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] robo_↔ FYIPProgram_FunctionPointerContainer)
- [CRobo_FYIPProgram_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetValve1](#) (unsigned char index, unsigned int valve1)
- unsigned int [GetValve1](#) (unsigned char index)
- void [SetValve2](#) (unsigned char index, unsigned int valve2)
- unsigned int [GetValve2](#) (unsigned char index)
- void [SetLength](#) (unsigned char index, int length)
- int [GetLength](#) (unsigned char index)
- void [Start](#) ()
- int [GetState](#) ()

Additional Inherited Members

11.98.1 Constructor & Destructor Documentation

11.98.1.1 CRobo_FYIPProgram_FunctionNet() [1/2] [CRobo_FYIPProgram_FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] robo_FYIPProgram_FunctionPointerContainer)

11.98.1.2 CRobo_FYIPProgram_FunctionNet() [2/2] [CRobo_FYIPProgram_FunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

11.98.2 Member Function Documentation

11.98.2.1 GetLength() `int GetLength (`
`unsigned char index)`

11.98.2.2 GetState() `int GetState ()`

11.98.2.3 GetValve1() `unsigned int GetValve1 (`
`unsigned char index)`

11.98.2.4 GetValve2() `unsigned int GetValve2 (`
`unsigned char index)`

11.98.2.5 SetLength() `void SetLength (`
`unsigned char index,`
`int length)`


11.98.2.6 SetValve1() `void SetValve1 (`
`unsigned char index,`
`unsigned int valve1)`

11.98.2.7 SetValve2() `void SetValve2 (`
`unsigned char index,`
`unsigned int valve2)`

11.98.2.8 Start() `void Start ()`

11.99 CRobo_FYITemp_FunctionNet Class Reference

Inheritance diagram for CRobo_FYITemp_FunctionNet:

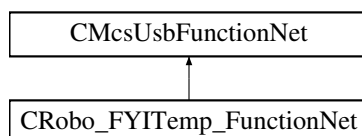

Public Member Functions

- [CRobo_FYITemp_FunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SetRegulatorOnOff](#) (unsigned char index, int onoff)
- int [GetRegulatorOnOff](#) (unsigned char index)
- void [SetSollTemp](#) (unsigned char index, int temp)
- int [GetSollTemp](#) (unsigned char index)
- void [SetPCoeff](#) (unsigned char index, int pcoeff)
- int [GetPCoeff](#) (unsigned char index)
- void [SetICoeff](#) (unsigned char index, int icoeff)
- int [GetICoeff](#) (unsigned char index)
- void [SetMaxPower](#) (unsigned char index, int power)
- int [GetMaxPower](#) (unsigned char index)

Additional Inherited Members

11.99.1 Constructor & Destructor Documentation

11.99.1.1 CRobo_FYITemp_FunctionNet() [CRobo_FYITemp_FunctionNet](#) (
 [CMcsUsbNet](#)[^] mcsusb)

11.99.2 Member Function Documentation

11.99.2.1 GetICoeff() int GetICoeff (
 unsigned char index)

11.99.2.2 GetMaxPower() int GetMaxPower (
 unsigned char index)

11.99.2.3 GetPCoeff() int GetPCoeff (
 unsigned char index)

11.99.2.4 GetRegulatorOnOff() int GetRegulatorOnOff (
 unsigned char index)

11.99.2.5 GetSollTemp() `int GetSollTemp (`
 `unsigned char index)`

11.99.2.6 SetICoeff() `void SetICoeff (`
 `unsigned char index,`
 `int icoeff)`

11.99.2.7 SetMaxPower() `void SetMaxPower (`
 `unsigned char index,`
 `int power)`


11.99.2.8 SetPCoeff() `void SetPCoeff (`
 `unsigned char index,`
 `int pcoeff)`

11.99.2.9 SetRegulatorOnOff() `void SetRegulatorOnOff (`
 `unsigned char index,`
 `int onoff)`

11.99.2.10 SetSollTemp() `void SetSollTemp (`
 `unsigned char index,`
 `int temp)`

11.100 CRoboDacqNet Class Reference

Inheritance diagram for CRoboDacqNet:

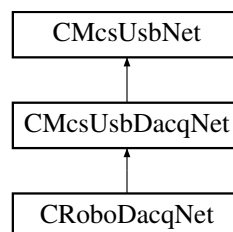

Public Member Functions

- [CRoboDacqNet](#) (void)
- [CRoboDacqNet](#) ([CRoboDeviceNet](#)^ robodevice)
- void [RunTable](#) ()
- void [RunTable](#) (int timeout)
- void [StopTable](#) ()
- void [StopTable](#) (int timeout)
- void [CancelTableLoop](#) ()
- void [CancelTableLoopAndStopTable](#) ()
- void [SetTriggerMaskValue](#) (unsigned int mask, unsigned int value, unsigned int virtualDevice)
- void [SetConfigurationBit](#) (unsigned short bit, bool value)
- void [SetConfigurationBitSupply](#) (bool value)
- void [SetConfigurationBitRelais](#) (bool value)
- void [SetConfigurationBitStream](#) (bool value)
- void [SetConfigurationBitAxc](#) (bool value)
- void [SetConfigurationBitCC_Gen](#) (bool value)
- void [SetConfigurationBitCV_Gen](#) (bool value)
- void [SetConfigurationBitRC_Gen](#) (bool value)
- void [SetConfigurationBitRV_Gen](#) (bool value)
- void [SetConfigurationBitBlu_Led](#) (bool value)
- void [SetConfigurationBitRed_Led](#) (bool value)
- void [SetConfigurationBitBlu_LedToggleSlow](#) (bool value)
- void [SetConfigurationBitRed_LedToggleSlow](#) (bool value)
- void [SetConfigurationBitBlu_LedToggleFast](#) (bool value)
- void [SetConfigurationBitRed_LedToggleFast](#) (bool value)
- void [SetConfigurationBitRed_LedSaturation](#) (bool value)
- void [SetSimulation](#) (unsigned int enable)
- void [SetUClamp](#) (int uClamp)
- void [SetIClamp](#) (int iClamp)
- void [SetPGain](#) (int pGain)
- void [SetIGain](#) (int iGain)
- void [SetFilter](#) (int filter)
- void [SetUVOffset](#) (int UVOffset)
- void [SetUCOffset](#) (int UCOffset)
- void [SetICOffset](#) (int ICOffset)
- void [SetCrossTalkOffset](#) (int CrossTalk)
- void [SetXGain](#) (int xGain)
- void [SetCrossTalkOptimum](#) (int cxOptimum)
- void [SetRecordingNumber](#) (unsigned int recordingNumber)
- void [ClampAmpRestart](#) ()
- void [DoRamp](#) (int startValue, int endValue, int duration, int mode)
- unsigned int [GetClampAmpSerialNumber](#) ()
- unsigned int [GetConfigurationBits](#) ()
- bool [GetConfigurationBit](#) (unsigned short bit)
- bool [GetConfigurationBitSupply](#) ()
- bool [GetConfigurationBitRelais](#) ()
- bool [GetConfigurationBitStream](#) ()
- bool [GetConfigurationBitAxc](#) ()
- bool [GetConfigurationBitCC_Gen](#) ()
- bool [GetConfigurationBitCV_Gen](#) ()
- bool [GetConfigurationBitRC_Gen](#) ()
- bool [GetConfigurationBitRV_Gen](#) ()
- bool [GetConfigurationBitBlu_Led](#) ()
- bool [GetConfigurationBitRed_Led](#) ()

- bool [GetConfigurationBitBlu_LedToggleSlow](#) ()
- bool [GetConfigurationBitRed_LedToggleSlow](#) ()
- bool [GetConfigurationBitBlu_LedToggleFast](#) ()
- bool [GetConfigurationBitRed_LedToggleFast](#) ()
- bool [GetConfigurationBitRed_LedSaturation](#) ()
- unsigned int [GetSimulation](#) ()
- int [GetUClamp](#) ()
- int [GetIClamp](#) ()
- int [GetPGain](#) ()
- int [GetIGain](#) ()
- int [GetFilter](#) ()
- int [GetUJOffset](#) ()
- int [GetUCOffset](#) ()
- int [GetICOffset](#) ()
- int [GetCrossTalkOffset](#) ()
- int [GetXGain](#) ()
- int [GetCrossTalkOptimum](#) ()
- unsigned int [GetRecordingNumber](#) ()
- int [GetResistanceC](#) ()
- int [GetResistanceV](#) ()
- int [GetCapacityC](#) ()
- int [GetCapacityV](#) ()
- int [GetCapacityX](#) ()
- int [GetUV](#) ()
- int [GetUC](#) ()
- int [GetIC](#) ()
- int [GetNUV_MS](#) ()
- int [GetNUC_MS](#) ()
- int [GetNIC_MS](#) ()
- void [SetAllDigout](#) (uint32_t value)
- uint32_t [GetAllDigout](#) ()
- void [SetCommand](#) (unsigned char command, int value)
- int [GetCommand](#) (unsigned char command)
- void [SetDigout](#) (uint16_t index, bool enable)
- bool [GetDigout](#) (uint16_t index)
- void [TableDefBegin](#) ()
- void [TableDefEnd](#) ()
- void [Table_Wait](#) (unsigned int tableWait)
- void [SetDownsampleFactor](#) (int index, int downsample_factor)
- void [SetFilterCoeffs](#) (int index, array< int > ^ coeffs)
- void [SetNoFilterCoeffs](#) (int index)
- int [GetDownsampleFactor](#) (int index)
- array< int > ^ [GetFilterCoeffs](#) (int index)
- void [Emu_SetElectrodeResists](#) (int emuElectrodeResist)
- void [Emu_SetCellResists](#) (int emuCellResist)
- void [Emu_SetCellCapacity](#) (int emuCellCapacity)
- void [Emu_SetCellPotential](#) (int emuCellPotential)
- void [Emu_SetNoise](#) (int emuNoise)
- int [Emu_GetElectrodeResists](#) ()
- int [Emu_GetCellResists](#) ()
- int [Emu_GetCellCapacity](#) ()
- int [Emu_GetCellPotential](#) ()
- int [Emu_GetNoise](#) ()
- void [SetDisplayText](#) (int index, String ^ displayText)
- void [SetScreen](#) (int screen)

- void [UpdateDisplay](#) ()
- String ^ [GetDisplayText](#) (int index)
- int [GetScreen](#) ()
- int [GetUpdateDisplay](#) ()

Static Public Attributes

- static const unsigned int [TriggerMask_Default](#) = 0xFF00
- static const unsigned int [TriggerValue_MoveAbs](#) = COMMAND_ROBO_MOVEABS
- static const unsigned int [TriggerValue_StartQueue](#) = COMMAND_ROBO_QUEUE
- static const unsigned int [VirtualDevice_ContinuousDacq](#) = 0
- static const unsigned int [VirtualDevice_TableRun](#) = 1

Additional Inherited Members

11.100.1 Constructor & Destructor Documentation

11.100.1.1 CRoboDacqNet() [1/2] [CRoboDacqNet](#) (
void)

11.100.1.2 CRoboDacqNet() [2/2] [CRoboDacqNet](#) (
[CRoboDeviceNet](#)^ *robodevice*)

11.100.2 Member Function Documentation

11.100.2.1 CancelTableLoop() void [CancelTableLoop](#) ()

11.100.2.2 CancelTableLoopAndStopTable() void [CancelTableLoopAndStopTable](#) ()

11.100.2.3 ClampAmpRestart() void [ClampAmpRestart](#) ()

11.100.2.4 DoRamp() `void DoRamp (`
 `int startValue,`
 `int endValue,`
 `int duration,`
 `int mode)`

11.100.2.5 Emu_GetCellCapacity() `int Emu_GetCellCapacity ()`

11.100.2.6 Emu_GetCellPotential() `int Emu_GetCellPotential ()`

11.100.2.7 Emu_GetCellResists() `int Emu_GetCellResists ()`

11.100.2.8 Emu_GetElectrodeResists() `int Emu_GetElectrodeResists ()`

11.100.2.9 Emu_GetNoise() `int Emu_GetNoise ()`

11.100.2.10 Emu_SetCellCapacity() `void Emu_SetCellCapacity (`
 `int emuCellCapacity)`

11.100.2.11 Emu_SetCellPotential() `void Emu_SetCellPotential (`
 `int emuCellPotential)`

11.100.2.12 Emu_SetCellResists() `void Emu_SetCellResists (`
 `int emuCellResist)`

11.100.2.13 Emu_SetElectrodeResists() `void Emu_SetElectrodeResists (`
 `int emuElectrodeResist)`

11.100.2.14 Emu_SetNoise() void Emu_SetNoise (
int *emuNoise*)

11.100.2.15 GetAllDigout() uint32_t GetAllDigout ()

11.100.2.16 GetCapacityC() int GetCapacityC ()

11.100.2.17 GetCapacityV() int GetCapacityV ()

11.100.2.18 GetCapacityX() int GetCapacityX ()

11.100.2.19 GetClampAmpSerialNumber() unsigned int GetClampAmpSerialNumber ()

11.100.2.20 GetCommand() int GetCommand (
unsigned char *command*)

11.100.2.21 GetConfigurationBit() bool GetConfigurationBit (
unsigned short *bit*)

11.100.2.22 GetConfigurationBitAxc() bool GetConfigurationBitAxc ()

11.100.2.23 GetConfigurationBitBlu_Led() bool GetConfigurationBitBlu_Led ()

11.100.2.24 GetConfigurationBitBlu_LedToggleFast() bool GetConfigurationBitBlu_LedToggleFast ()

11.100.2.25 GetConfigurationBitBlu_LedToggleSlow() `bool GetConfigurationBitBlu_LedToggleSlow ()`

11.100.2.26 GetConfigurationBitCC_Gen() `bool GetConfigurationBitCC_Gen ()`

11.100.2.27 GetConfigurationBitCV_Gen() `bool GetConfigurationBitCV_Gen ()`

11.100.2.28 GetConfigurationBitRC_Gen() `bool GetConfigurationBitRC_Gen ()`

11.100.2.29 GetConfigurationBitRed_Led() `bool GetConfigurationBitRed_Led ()`

11.100.2.30 GetConfigurationBitRed_LedSaturation() `bool GetConfigurationBitRed_LedSaturation ()`

11.100.2.31 GetConfigurationBitRed_LedToggleFast() `bool GetConfigurationBitRed_LedToggleFast ()`

11.100.2.32 GetConfigurationBitRed_LedToggleSlow() `bool GetConfigurationBitRed_LedToggleSlow ()`

11.100.2.33 GetConfigurationBitRelais() `bool GetConfigurationBitRelais ()`

11.100.2.34 GetConfigurationBitRV_Gen() `bool GetConfigurationBitRV_Gen ()`

11.100.2.35 GetConfigurationBits() `unsigned int GetConfigurationBits ()`

11.100.2.36 GetConfigurationBitStream() `bool GetConfigurationBitStream ()`

11.100.2.37 GetConfigurationBitSupply() `bool GetConfigurationBitSupply ()`

11.100.2.38 GetCrossTalkOffset() `int GetCrossTalkOffset ()`

11.100.2.39 GetCrossTalkOptimum() `int GetCrossTalkOptimum ()`

11.100.2.40 GetDigout() `bool GetDigout (`
`uint16_t index)`

11.100.2.41 GetDisplayText() `String ^ GetDisplayText (`
`int index)`

11.100.2.42 GetDownsampleFactor() `int GetDownsampleFactor (`
`int index)`

11.100.2.43 GetFilter() `int GetFilter ()`

11.100.2.44 GetFilterCoeffs() `array<int> ^ GetFilterCoeffs (`
`int index)`

11.100.2.45 GetIC() `int GetIC ()`

11.100.2.46 GetIClamp() `int GetIClamp ()`

11.100.2.47 GetICOffset() `int GetICOffset ()`

11.100.2.48 GetIGain() `int GetIGain ()`

11.100.2.49 GetNIC_MS() `int GetNIC_MS ()`

11.100.2.50 GetNUC_MS() `int GetNUC_MS ()`

11.100.2.51 GetNUV_MS() `int GetNUV_MS ()`

11.100.2.52 GetPGain() `int GetPGain ()`

11.100.2.53 GetRecordingNumber() `unsigned int GetRecordingNumber ()`

11.100.2.54 GetResistanceC() `int GetResistanceC ()`

11.100.2.55 GetResistanceV() `int GetResistanceV ()`

11.100.2.56 GetScreen() `int GetScreen ()`

11.100.2.57 GetSimulation() `unsigned int GetSimulation ()`

11.100.2.58 GetUC() `int GetUC ()`

11.100.2.59 GetUClamp() `int GetUClamp ()`

11.100.2.60 GetUCOffset() `int GetUCOffset ()`

11.100.2.61 GetUpdateDisplay() `int GetUpdateDisplay ()`

11.100.2.62 GetUV() `int GetUV ()`

11.100.2.63 GetUVOffset() `int GetUVOffset ()`

11.100.2.64 GetXGain() `int GetXGain ()`

11.100.2.65 RunTable() [1/2] `void RunTable ()`

11.100.2.66 RunTable() [2/2] `void RunTable (`
`int timeout)`

11.100.2.67 SetAllDigout() `void SetAllDigout (`
`uint32_t value)`

11.100.2.68 SetCommand() `void SetCommand (`
`unsigned char command,`
`int value)`

11.100.2.69 SetConfigurationBit() void SetConfigurationBit (
 unsigned short *bit*,
 bool *value*)

11.100.2.70 SetConfigurationBitAxc() void SetConfigurationBitAxc (
 bool *value*)

11.100.2.71 SetConfigurationBitBlu_Led() void SetConfigurationBitBlu_Led (
 bool *value*)

11.100.2.72 SetConfigurationBitBlu_LedToggleFast() void SetConfigurationBitBlu_LedToggleFast (
 bool *value*)

11.100.2.73 SetConfigurationBitBlu_LedToggleSlow() void SetConfigurationBitBlu_LedToggleSlow (
 bool *value*)

11.100.2.74 SetConfigurationBitCC_Gen() void SetConfigurationBitCC_Gen (
 bool *value*)

11.100.2.75 SetConfigurationBitCV_Gen() void SetConfigurationBitCV_Gen (
 bool *value*)

11.100.2.76 SetConfigurationBitRC_Gen() void SetConfigurationBitRC_Gen (
 bool *value*)

11.100.2.77 SetConfigurationBitRed_Led() void SetConfigurationBitRed_Led (
 bool *value*)

11.100.2.78 SetConfigurationBitRed_LedSaturation() void SetConfigurationBitRed_LedSaturation (bool *value*)

11.100.2.79 SetConfigurationBitRed_LedToggleFast() void SetConfigurationBitRed_LedToggleFast (bool *value*)

11.100.2.80 SetConfigurationBitRed_LedToggleSlow() void SetConfigurationBitRed_LedToggleSlow (bool *value*)

11.100.2.81 SetConfigurationBitRelais() void SetConfigurationBitRelais (bool *value*)

11.100.2.82 SetConfigurationBitRV_Gen() void SetConfigurationBitRV_Gen (bool *value*)

11.100.2.83 SetConfigurationBitStream() void SetConfigurationBitStream (bool *value*)

11.100.2.84 SetConfigurationBitSupply() void SetConfigurationBitSupply (bool *value*)

11.100.2.85 SetCrossTalkOffset() void SetCrossTalkOffset (int *CrossTalk*)

11.100.2.86 SetCrossTalkOptimum() void SetCrossTalkOptimum (int *cxOptimum*)

11.100.2.87 SetDigout() void SetDigout (
 uint16_t *index*,
 bool *enable*)

11.100.2.88 SetDisplayText() void SetDisplayText (
 int *index*,
 String^ *displayText*)

11.100.2.89 SetDownsampleFactor() void SetDownsampleFactor (
 int *index*,
 int *downsample_factor*)

11.100.2.90 SetFilter() void SetFilter (
 int *filter*)

11.100.2.91 SetFilterCoeffs() void SetFilterCoeffs (
 int *index*,
 array< int >^ *coeffs*)

11.100.2.92 SetIClamp() void SetIClamp (
 int *iClamp*)

11.100.2.93 SetICOffset() void SetICOffset (
 int *ICOffset*)

11.100.2.94 SetIGain() void SetIGain (
 int *iGain*)

11.100.2.95 SetNoFilterCoeffs() void SetNoFilterCoeffs (
 int *index*)

11.100.2.96 SetPGain() void SetPGain (
int *pGain*)

11.100.2.97 SetRecordingNumber() void SetRecordingNumber (
unsigned int *recordingNumber*)

11.100.2.98 SetScreen() void SetScreen (
int *screen*)

11.100.2.99 SetSimulation() void SetSimulation (
unsigned int *enable*)

11.100.2.100 SetTriggerMaskValue() void SetTriggerMaskValue (
unsigned int *mask*,
unsigned int *value*,
unsigned int *virtualDevice*)

11.100.2.101 SetUClamp() void SetUClamp (
int *uClamp*)

11.100.2.102 SetUCOffset() void SetUCOffset (
int *UCOffset*)

11.100.2.103 SetUVOffset() void SetUVOffset (
int *UVOffset*)

11.100.2.104 SetXGain() void SetXGain (
int *xGain*)

11.100.2.105 StopTable() [1/2] `void StopTable ()`

11.100.2.106 StopTable() [2/2] `void StopTable (`
`int timeout)`

11.100.2.107 Table_Wait() `void Table_Wait (`
`unsigned int tableWait)`

11.100.2.108 TableDefBegin() `void TableDefBegin ()`

11.100.2.109 TableDefEnd() `void TableDefEnd ()`

11.100.2.110 UpdateDisplay() `void UpdateDisplay ()`

11.100.3 Member Data Documentation

11.100.3.1 TriggerMask_Default `const unsigned int TriggerMask_Default = 0xFF00 [static]`

11.100.3.2 TriggerValue_MoveAbs `const unsigned int TriggerValue_MoveAbs = COMMAND_ROBO_↔`
`MOVEABS [static]`

11.100.3.3 TriggerValue_StartQueue `const unsigned int TriggerValue_StartQueue = COMMAND_ROBO_↔`
`QUEUE [static]`


11.100.3.4 VirtualDevice_ContinousDacq `const unsigned int VirtualDevice_ContinousDacq = 0`
`[static]`

11.100.3.5 VirtualDevice_TableRun `const unsigned int VirtualDevice_TableRun = 1 [static]`

11.101 CRoboDeviceNet Class Reference

[CRoboDeviceNet](#) is the base class for all Robo platform based devices

Inheritance diagram for CRoboDeviceNet:

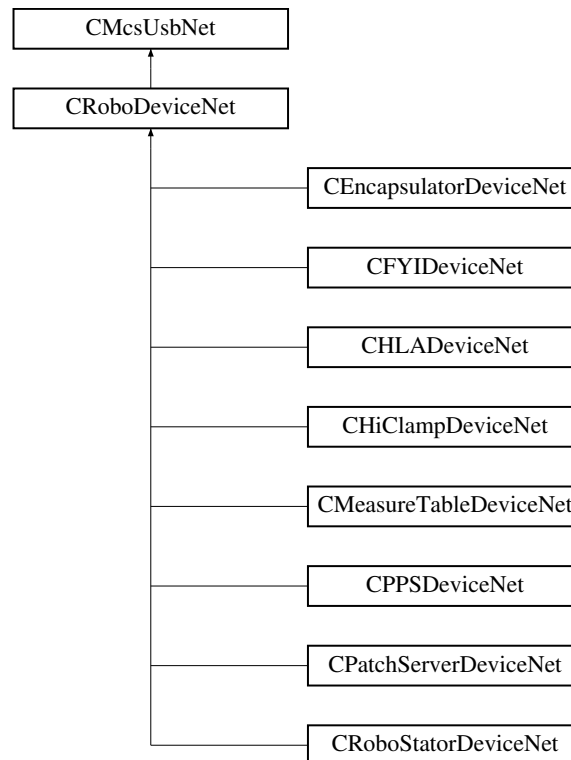

Classes

- class [RoboMainLowLevelCommands](#)

Public Member Functions

- [CRoboDeviceNet](#) (void)
- [~CRoboDeviceNet](#) (void)
- void [SetInMovement](#) ()
Low level command, sets the internal state to "In Movement"
- bool [GetInMovement](#) ()
Low level command, gets the internal state "In Movement"
- uint32_t [GetMovementError](#) ()
Low level command, gets the error of the last movement end
- void [FindReference](#) (unsigned char busaddress, char axes)
- void [FindReference](#) (unsigned char busaddress, char axes, int timeout)
Searches the reference position of the motor
- void [MoveAbs](#) (unsigned char busaddress, char axes, int x, int y)
Moves the motor to the new absolute position

- void [MoveAbs](#) (unsigned char busaddress, char axes, int x, int y, int timeout)
Moves the motor to the new absolute position
- void [MoveAbs](#) (unsigned char busaddress, char axes, array< int >^ pos)
Moves the motor to the new absolute position
- void [MoveAbs](#) (unsigned char busaddress, char axes, array< int >^ pos, int timeout)
Moves the motor to the new absolute position
- void [StopMovement](#) (unsigned char busaddress, char axes)
- void [StopMovement](#) (unsigned char busaddress, char axes, int timeout)
Stops the current movement
- void [SetCurrentAndAir](#) (unsigned char busaddress, char axes, unsigned short onoff)
- void [SetCurrentAndAir](#) (unsigned char busaddress, char axes, unsigned short onoff, int timeout)
- bool [IsQueueEnabled](#) ()
- void [EnableQueue](#) (bool enable)
- bool [IsQueueStarted](#) ()
- void [StartQueue](#) (bool start)
- void [WaitTimer](#) (uint32_t waittime, int timeout)
- void [CancelPoolLoop](#) ()
- void [CancelPoolLoopAndStopMovement](#) ()
- void [GetCurrentPosition](#) (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out]int% x, [System::Runtime::InteropServices::Out]int% y)
Gets the current position of motors
- void [SetAirValve](#) (unsigned int onoff)
- unsigned int [GetAirValve](#) ()
- void [NullCommand](#) (unsigned int marker)
- unsigned int [GetVoltageValves](#) ()
- unsigned int [GetVoltageRs485A](#) ()
- unsigned int [GetVoltageRs485B](#) ()
- unsigned int [GetVoltageAirvalve](#) ()
- unsigned int [GetCurrentAirvalve](#) ()
- unsigned int [GetVoltage12V](#) ()
- unsigned int [GetAirpressure](#) ()
- unsigned int [GetVoltage5V](#) ()
- unsigned int [GetErrorVoltageValves](#) ()
- unsigned int [GetErrorVoltageRs485A](#) ()
- unsigned int [GetErrorVoltageRs485B](#) ()
- unsigned int [GetErrorVoltageAirvalve](#) ()
- unsigned int [GetErrorCurrentAirvalve](#) ()
- unsigned int [GetErrorVoltage12V](#) ()
- unsigned int [GetErrorAirpressure](#) ()
- unsigned int [GetErrorVoltage5V](#) ()
- void [SetVoltageValvesLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageRs485ALimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageRs485BLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetVoltageAirvalveLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetCurrentAirvalveLimit](#) (unsigned int lowercurrent, unsigned int uppercurrent)
- void [SetVoltage12VLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [SetAirpressureLimit](#) (unsigned int lowerpressure, unsigned int upperpressure)
- void [SetVoltage5VLimit](#) (unsigned int lowervoltage, unsigned int uppervoltage)
- void [GetVoltageValvesLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetVoltageRs485ALimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)
- void [GetVoltageRs485BLimit](#) ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System::Runtime::InteropServices::Out] unsigned int% uppervoltage)

- void [GetVoltageAirvalveLimit](#) ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void [GetCurrentAirvalveLimit](#) ([System::Runtime::InteropServices::Out]unsigned int% lowercurrent, [System::Runtime::InteropServices::Out]unsigned int% uppercurrent)
- void [GetVoltage12VLimit](#) ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void [GetAirpressureLimit](#) ([System::Runtime::InteropServices::Out]unsigned int% lowerpressure, [System::Runtime::InteropServices::Out]unsigned int% upperpressure)
- void [GetVoltage5VLimit](#) ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void [SetMinPressure](#) (int pressure)
- int [GetMinPressure](#) ()

Static Public Attributes

- static const uint32_t [RoboError_Base](#) = (0xA0110000L)
- static const uint32_t [RoboError_UnknownCommand](#) = ((0xA0110000L))
- static const uint32_t [RoboError_Timeout](#) = ((0xA0110000L) | 0x0001)
- static const uint32_t [RoboError_Pressure](#) = ((0xA0110000L) | 0x0002)
- static const uint32_t [RoboError_RangeExceeded](#) = ((0xA0110000L) | 0x0003)
- static const uint32_t [RoboError_CommunicationTimeout](#) = ((0xA0110000L) | 0x0004)
- static const uint32_t [RoboError_AnotherMaster](#) = ((0xA0110000L) | 0x0005)
- static const uint32_t [RoboError_FindReferenceMethod](#) = ((0xA0110000L) | 0x0006)
- static const uint32_t [RoboError_NoSpeedOrAcceleration](#) = ((0xA0110000L) | 0x0007)
- static const uint32_t [RoboError_NoEndSwitch](#) = ((0xA0110000L) | 0x0008)
- static const uint32_t [RoboError_CannotEscapeEndSwitch](#) = ((0xA0110000L) | 0x0009)
- static const uint32_t [RoboError_CommandAlreadyInProgress](#) = ((0xA0110000L) | 0x000A)
- static const uint32_t [RoboError_NoReference](#) = ((0xA0110000L) | 0x000B)
- static const uint32_t [RoboError_OverPressure](#) = ((0xA0110000L) | 0x000C)
- static const uint32_t [RoboError_Phase0OutOfRange](#) = ((0xA0110000L) | 0x000D)
- static const uint32_t [RoboError_PeristalticTimeout](#) = ((0xA0110000L) | 0x000E)
- static const uint32_t [RoboError_GilsonTimeout](#) = ((0xA0110000L) | 0x000F)
- static const uint32_t [RoboError_GilsonWrondID](#) = ((0xA0110000L) | 0x0010)
- static const uint32_t [RoboError_GilsonCommandPending](#) = ((0xA0110000L) | 0x0011)
- static const uint32_t [RoboError_ParameterNotAllowed](#) = ((0xA0110000L) | 0x0012)
- static const uint32_t [RoboError_StateChangeNotPossible](#) = ((0xA0110000L) | 0x0013)
- static const uint32_t [RoboError_CommandNotPossible](#) = ((0xA0110000L) | 0x0014)
- static const uint32_t [RoboError_DacqNotReady](#) = ((0xA0110000L) | 0x0015)
- static const uint32_t [RoboError_NoMoreData](#) = ((0xA0110000L) | 0x0016)
- static const uint32_t [RoboError_McsBus_UnknownCommand](#) = ((0xA0110000L) | 0x003F)
- static const uint32_t [RoboError_DLLMovementTimeout](#) = ((0xA0110000L) | 0x1001)
- static const uint32_t [RoboError_PollLoopCanceled](#) = ((0xA0110000L) | 0x1002)
- static const uint32_t [RoboError_PollLoopCanceledAndStopMovement](#) = ((0xA0110000L) | 0x1003)
- static const byte [McsBus_XY](#) = 1
McsBus address for the xy-plane
- static const byte [McsBus_ZI](#) = 2
McsBus address for the z and i axes
- static const byte [Axis_X](#) = 0
Axis number of x for axis argument
- static const byte [Axis_Y](#) = 1
Axis number of y for axis argument
- static const byte [Axis_Z](#) = 0
Axis number of z for axis argument

- static const byte [Axis_I](#) = 1
Axis number of i for axis argument
- static const char [Axes_X](#) = 1
Bit pattern for x axis for axes argument
- static const char [Axes_Y](#) = 2
Bit pattern for y axis for axes argument
- static const char [Axes_Z](#) = 1
Bit pattern for z axis for axes argument
- static const char [Axes_I](#) = 2
Bit pattern for i axis for axes argument

Properties

- [CMcsBusNet](#)[^] [McsBus](#) [get]
- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]
- [RoboMainLowLevelCommands](#)[^] [RoboMainLowLevelCommand](#) [get]

Events

- [RoboStatusEventDelegate](#)[^] [RoboStatusEvent](#)

Additional Inherited Members

11.101.1 Detailed Description

[CRoboDeviceNet](#) is the base class for all Robo platform based devices

11.101.2 Constructor & Destructor Documentation

11.101.2.1 [CRoboDeviceNet](#)() [CRoboDeviceNet](#) (
void)

11.101.2.2 [~CRoboDeviceNet](#)() [~CRoboDeviceNet](#) (
void)

11.101.3 Member Function Documentation

11.101.3.1 [CancelPoolLoop](#)() void [CancelPoolLoop](#) ()

11.101.3.2 CancelPoolLoopAndStopMovement() `void CancelPoolLoopAndStopMovement ()`

11.101.3.3 EnableQueue() `void EnableQueue (`
`bool enable)`

11.101.3.4 FindReference() [1/2] `void FindReference (`
`unsigned char busaddress,`
`char axes)`

11.101.3.5 FindReference() [2/2] `void FindReference (`
`unsigned char busaddress,`
`char axes,`
`int timeout)`

Searches the reference position of the motor

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

11.101.3.6 GetAirpressure() `unsigned int GetAirpressure ()`

11.101.3.7 GetAirpressureLimit() `void GetAirpressureLimit (`
`[System::Runtime::InteropServices::Out] unsigned int% lowerpressure,`
`[System::Runtime::InteropServices::Out] unsigned int% upperpressure)`

11.101.3.8 GetAirValve() `unsigned int GetAirValve ()`

11.101.3.9 GetCurrentAirvalve() `unsigned int GetCurrentAirvalve ()`

11.101.3.10 GetCurrentAirvalveLimit() void GetCurrentAirvalveLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowercurrent*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppercurrent*)

11.101.3.11 GetCurrentPosition() void GetCurrentPosition (
 unsigned char *busaddress*,
 char *axes*,
 [System::Runtime::InteropServices::Out] int% *x*,
 [System::Runtime::InteropServices::Out] int% *y*)

Gets the current position of motors

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>x</i>	Current position of first axis if pattern in axes is set
<i>y</i>	Current position of second axis if pattern in axes is set

11.101.3.12 GetErrorAirpressure() unsigned int GetErrorAirpressure ()

11.101.3.13 GetErrorCurrentAirvalve() unsigned int GetErrorCurrentAirvalve ()

11.101.3.14 GetErrorVoltage12V() unsigned int GetErrorVoltage12V ()

11.101.3.15 GetErrorVoltage5V() unsigned int GetErrorVoltage5V ()

11.101.3.16 GetErrorVoltageAirvalve() unsigned int GetErrorVoltageAirvalve ()

11.101.3.17 GetErrorVoltageRs485A() unsigned int GetErrorVoltageRs485A ()

11.101.3.18 GetErrorVoltageRs485B() unsigned int GetErrorVoltageRs485B ()

11.101.3.19 GetErrorVoltageValves() unsigned int GetErrorVoltageValves ()

11.101.3.20 GetInMovement() bool GetInMovement ()

Low level command, gets the internal state "In Movement"

11.101.3.21 GetMinPressure() int GetMinPressure ()

11.101.3.22 GetMovementError() uint32_t GetMovementError ()

Low level command, gets the error of the last movement end

11.101.3.23 GetVoltage12V() unsigned int GetVoltage12V ()

11.101.3.24 GetVoltage12VLimit() void GetVoltage12VLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.25 GetVoltage5V() unsigned int GetVoltage5V ()

11.101.3.26 GetVoltage5VLimit() void GetVoltage5VLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.27 GetVoltageAirvalve() unsigned int GetVoltageAirvalve ()

11.101.3.28 GetVoltageAirvalveLimit() void GetVoltageAirvalveLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.29 GetVoltageRs485A() unsigned int GetVoltageRs485A ()

11.101.3.30 GetVoltageRs485ALimit() void GetVoltageRs485ALimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.31 GetVoltageRs485B() unsigned int GetVoltageRs485B ()

11.101.3.32 GetVoltageRs485BLimit() void GetVoltageRs485BLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.33 GetVoltageValves() unsigned int GetVoltageValves ()

11.101.3.34 GetVoltageValvesLimit() void GetVoltageValvesLimit (
 [System::Runtime::InteropServices::Out] unsigned int% *lowervoltage*,
 [System::Runtime::InteropServices::Out] unsigned int% *uppervoltage*)

11.101.3.35 IsQueueEnabled() bool IsQueueEnabled ()

11.101.3.36 IsQueueStarted() bool IsQueueStarted ()

11.101.3.37 MoveAbs() [1/4] void MoveAbs (
 unsigned char *busaddress*,
 char *axes*,
 array< int >^ *pos*)

Moves the motor to the new absolute position

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>pos</i>	Positions of the axis 0 to 3, if pattern in axes is set

11.101.3.38 MoveAbs() [2/4] `void MoveAbs (`
 unsigned char *busaddress*,
 char *axes*,
 array< int >^ *pos*,
 int *timeout*)

Moves the motor to the new absolute position

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>pos</i>	Positions of the axis 0 to 3, if pattern in axes is set
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

11.101.3.39 MoveAbs() [3/4] `void MoveAbs (`
 unsigned char *busaddress*,
 char *axes*,
 int *x*,
 int *y*)

Moves the motor to the new absolute position

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>x</i>	Position of first axis, if pattern in axes is set
<i>y</i>	Position of second axis if pattern in axes is set

11.101.3.40 MoveAbs() [4/4] `void MoveAbs (`
 unsigned char *busaddress*,
 char *axes*,
 int *x*,
 int *y*,
 int *timeout*)

Moves the motor to the new absolute position

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>x</i>	Position of first axis, if pattern in axes is set
<i>y</i>	Position of second axis if pattern in axes is set
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

11.101.3.41 NullCommand() `void NullCommand (`
 `unsigned int marker)`

11.101.3.42 SetAirpressureLimit() `void SetAirpressureLimit (`
 `unsigned int lowerpressure,`
 `unsigned int upperpressure)`

11.101.3.43 SetAirValve() `void SetAirValve (`
 `unsigned int onoff)`

11.101.3.44 SetCurrentAirvalveLimit() `void SetCurrentAirvalveLimit (`
 `unsigned int lowercurrent,`
 `unsigned int uppercurrent)`

11.101.3.45 SetCurrentAndAir() [1/2] `void SetCurrentAndAir (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned short onoff)`

11.101.3.46 SetCurrentAndAir() [2/2] `void SetCurrentAndAir (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned short onoff,`
 `int timeout)`

11.101.3.47 SetInMovement() `void SetInMovement ()`

Low level command, sets the internal state to "In Movement"

11.101.3.48 SetMinPressure() `void SetMinPressure (`
`int pressure)`

11.101.3.49 SetVoltage12VLimit() `void SetVoltage12VLimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.50 SetVoltage5VLimit() `void SetVoltage5VLimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.51 SetVoltageAirvalveLimit() `void SetVoltageAirvalveLimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.52 SetVoltageRs485ALimit() `void SetVoltageRs485ALimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.53 SetVoltageRs485BLimit() `void SetVoltageRs485BLimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.54 SetVoltageValvesLimit() `void SetVoltageValvesLimit (`
`unsigned int lowervoltage,`
`unsigned int uppervoltage)`

11.101.3.55 StartQueue() `void StartQueue (`
 `bool start)`

11.101.3.56 StopMovement() [1/2] `void StopMovement (`
 `unsigned char busaddress,`
 `char axes)`

11.101.3.57 StopMovement() [2/2] `void StopMovement (`
 `unsigned char busaddress,`
 `char axes,`
 `int timeout)`

Stops the current movement

Parameters

<i>busaddress</i>	Address of the McsBus
<i>axes</i>	Bit pattern of axes to drive
<i>timeout</i>	Timeout of maximal waiting for the end of the command (-1 is forever)

11.101.3.58 WaitTimer() `void WaitTimer (`
 `uint32_t waittime,`
 `int timeout)`

11.101.4 Member Data Documentation

11.101.4.1 Axes_I `const char Axes_I = 2 [static]`

Bit pattern for i axis for axes argument

11.101.4.2 Axes_X `const char Axes_X = 1 [static]`

Bit pattern for x axis for axes argument

11.101.4.3 Axes_Y `const char Axes_Y = 2 [static]`

Bit pattern for y axis for axes argument

11.101.4.4 Axes_Z `const char Axes_Z = 1 [static]`

Bit pattern for z axis for axes argument

11.101.4.5 Axis_I `const byte Axis_I = 1 [static]`

Axis number of i for axis argument

11.101.4.6 Axis_X `const byte Axis_X = 0 [static]`

Axis number of x for axis argument

11.101.4.7 Axis_Y `const byte Axis_Y = 1 [static]`

Axis number of y for axis argument

11.101.4.8 Axis_Z `const byte Axis_Z = 0 [static]`

Axis number of z for axis argument

11.101.4.9 McsBus_XY `const byte McsBus_XY = 1 [static]`

McsBus address for the xy-plane

11.101.4.10 McsBus_ZI `const byte McsBus_ZI = 2 [static]`

McsBus address for the z and i axes

11.101.4.11 RoboError_AnotherMaster `const uint32_t RoboError_AnotherMaster = ((0xA0110000L) | 0x0005) [static]`

11.101.4.12 RoboError_Base `const uint32_t RoboError_Base = (0xA0110000L) [static]`

11.101.4.13 RoboError_CannotEscapeEndSwitch `const uint32_t RoboError_CannotEscapeEndSwitch = ((0xA0110000L) | 0x0009) [static]`

11.101.4.14 RoboError_CommandAlreadyInProgress `const uint32_t RoboError_CommandAlreadyInProgress = ((0xA0110000L) | 0x000A) [static]`

11.101.4.15 RoboError_CommandNotPossible `const uint32_t RoboError_CommandNotPossible = ((0xA0110000L) | 0x0014) [static]`

11.101.4.16 RoboError_CommunicationTimeout `const uint32_t RoboError_CommunicationTimeout = ((0xA0110000L) | 0x0004) [static]`

11.101.4.17 RoboError_DacqNotReady `const uint32_t RoboError_DacqNotReady = ((0xA0110000L) | 0x0015) [static]`

11.101.4.18 RoboError_DLLMovementTimeout `const uint32_t RoboError_DLLMovementTimeout = ((0xA0110000L) | 0x1001) [static]`

11.101.4.19 RoboError_FindReferenceMethod `const uint32_t RoboError_FindReferenceMethod = ((0xA0110000L) | 0x0006) [static]`

11.101.4.20 RoboError_GilsonCommandPending `const uint32_t RoboError_GilsonCommandPending = ((0xA0110000L) | 0x0011) [static]`

11.101.4.21 RoboError_GilsonTimeout `const uint32_t RoboError_GilsonTimeout = ((0xA0110000L) | 0x000F) [static]`

11.101.4.22 RoboError_GilsonWrondID `const uint32_t RoboError_GilsonWrondID = ((0xA0110000L) | 0x0010) [static]`

11.101.4.23 RoboError_McsBus_UnknownCommand `const uint32_t RoboError_McsBus_UnknownCommand = ((0xA0110000L) | 0x003F) [static]`

11.101.4.24 RoboError_NoEndSwitch `const uint32_t RoboError_NoEndSwitch = ((0xA0110000L) | 0x0008) [static]`

11.101.4.25 RoboError_NoMoreData `const uint32_t RoboError_NoMoreData = ((0xA0110000L) | 0x0016) [static]`

11.101.4.26 RoboError_NoReference `const uint32_t RoboError_NoReference = ((0xA0110000L) | 0x000B) [static]`

11.101.4.27 RoboError_NoSpeedOrAcceleration `const uint32_t RoboError_NoSpeedOrAcceleration = ((0xA0110000L) | 0x0007) [static]`

11.101.4.28 RoboError_OverPressure `const uint32_t RoboError_OverPressure = ((0xA0110000L) | 0x000C) [static]`

11.101.4.29 RoboError_ParameterNotAllowed `const uint32_t RoboError_ParameterNotAllowed = ((0xA0110000L) | 0x0012) [static]`

11.101.4.30 RoboError_PeristalticTimeout `const uint32_t RoboError_PeristalticTimeout = ((0xA0110000L) | 0x000E) [static]`

11.101.4.31 RoboError_Phase0OutOfRange `const uint32_t RoboError_Phase0OutOfRange = ((0x↵A0110000L) | 0x000D) [static]`

11.101.4.32 RoboError_PollLoopCanceled `const uint32_t RoboError_PollLoopCanceled = ((0x↵A0110000L) | 0x1002) [static]`

11.101.4.33 RoboError_PollLoopCanceledAndStopMovement `const uint32_t RoboError_PollLoop↵CanceledAndStopMovement = ((0xA0110000L) | 0x1003) [static]`

11.101.4.34 RoboError_Pressure `const uint32_t RoboError_Pressure = ((0xA0110000L) | 0x0002) [static]`

11.101.4.35 RoboError_RangeExceeded `const uint32_t RoboError_RangeExceeded = ((0xA0110000L) | 0x0003) [static]`

11.101.4.36 RoboError_StateChangeNotPossible `const uint32_t RoboError_StateChangeNotPossible = ((0xA0110000L) | 0x0013) [static]`

11.101.4.37 RoboError_Timeout `const uint32_t RoboError_Timeout = ((0xA0110000L) | 0x0001) [static]`

11.101.4.38 RoboError_UnknownCommand `const uint32_t RoboError_UnknownCommand = ((0x↵A0110000L)) [static]`

11.101.5 Property Documentation

11.101.5.1 McsBus `CMcsBusNet^ McsBus [get]`

11.101.5.2 McsBus_MotorControl [CMcsBus_MotorControlNet](#)[^] McsBus_MotorControl [get]


11.101.5.3 RoboMainLowLevelCommand [RoboMainLowLevelCommands](#)[^] RoboMainLowLevelCommand [get]

11.101.6 Event Documentation

11.101.6.1 RoboStatusEvent [RoboStatusEventDelegate](#)[^] RoboStatusEvent

11.102 CRoboFluidDeviceNet Class Reference

Inheritance diagram for CRoboFluidDeviceNet:

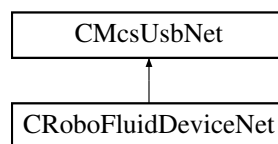

Public Member Functions

- [CRoboFluidDeviceNet](#) (void)
- [~CRoboFluidDeviceNet](#) (void)
- void [SetValve](#) (int value)
Open or Close valves.
- void [SetSingleValve](#) (int valve, bool onoff)
Opens or Closes a valve.
- int [GetValve](#) ()
Query the state of the values.
- bool [GetSingleValve](#) (int valve)
Query the state of a valve.
- void [CloseAllValves](#) ()
- void [PumpOn](#) (int index, short speed)
- void [SetPumpSpeed](#) (int index, short speed)
- void [PumpOff](#) (int index)
- short [GetPumpSpeed](#) (int index)
- bool [IsPumpMotorOn](#) (int index)

Protected Attributes

- CRoboFluidDevice * [m_pRoboFluidDevice](#)
- [CMcsBus_MotorControlNet](#)[^] [m_pMcsBus_MotorControlNet](#)

Properties

- [CMcsBus_MotorControlNet](#)[^] [McsBus_MotorControl](#) [get]

Additional Inherited Members

11.102.1 Constructor & Destructor Documentation

11.102.1.1 CRoboFluidDeviceNet() [CRoboFluidDeviceNet](#) (
 void)

11.102.1.2 ~CRoboFluidDeviceNet() [~CRoboFluidDeviceNet](#) (
 void)

11.102.2 Member Function Documentation

11.102.2.1 CloseAllValves() void CloseAllValves ()

11.102.2.2 GetPumpSpeed() short GetPumpSpeed (
 int *index*)

11.102.2.3 GetSingleValve() bool GetSingleValve (
 int *valve*)

Query the state of a valve.

Parameters

<i>valve</i>	number of valve /*!
--------------	---------------------

Returns

state of the valve

11.102.2.4 GetValve() `int GetValve ()`

Query the state of the values.

Returns

the current state of the valves as a bit pattern.

11.102.2.5 IsPumpMotorOn() `bool IsPumpMotorOn (`
`int index)`**11.102.2.6 PumpOff()** `void PumpOff (`
`int index)`**11.102.2.7 PumpOn()** `void PumpOn (`
`int index,`
`short speed)`**11.102.2.8 SetPumpSpeed()** `void SetPumpSpeed (`
`int index,`
`short speed)`**11.102.2.9 SetSingleValve()** `void SetSingleValve (`
`int valve,`
`bool onoff)`

Opens or Closes a valve.

Parameters

<i>valve</i>	number of valve to be changed /*!
<i>onoff</i>	open or close the valve

11.102.2.10 SetValve() `void SetValve (`
`int value)`

Open or Close valves.

Parameters

<i>value</i>	bit pattern of valves which should be open.
--------------	---

11.102.3 Member Data Documentation

11.102.3.1 m_pMcsBus_MotorControlNet [CMcsBus_MotorControlNet](#) ^ m_pMcsBus_MotorControlNet
[protected]

11.102.3.2 m_pRoboFluidDevice [CRoboFluidDevice*](#) m_pRoboFluidDevice [protected]


11.102.4 Property Documentation

11.102.4.1 McsBus_MotorControl [CMcsBus_MotorControlNet](#) ^ McsBus_MotorControl [get]

11.103 CRoboInjectDeviceNet Class Reference

[CRoboInjectDeviceNet](#) is the to control the MCS RoboInject device

Inheritance diagram for CRoboInjectDeviceNet:

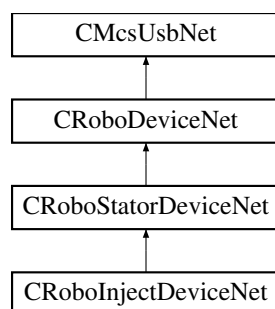

Public Member Functions

- [CRoboInjectDeviceNet](#) (void)

Additional Inherited Members

11.103.1 Detailed Description

[CRobolInjectDeviceNet](#) is the to control the MCS RobolInject device


11.103.2 Constructor & Destructor Documentation

11.103.2.1 CRobolInjectDeviceNet() [CRobolInjectDeviceNet](#) (
void)

11.104 CRoboocyte2DeviceNet Class Reference

[CRoboocyte2DeviceNet](#) is the class to control the MCS Roboocyte2 device

Inheritance diagram for CRoboocyte2DeviceNet:

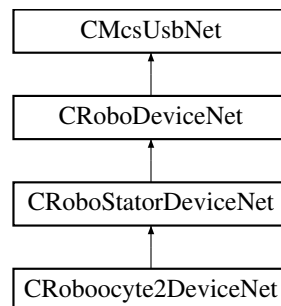

Public Member Functions

- [CRoboocyte2DeviceNet](#) (void)
- void [SetAxisLED](#) (bool onoff)
- bool [GetAxisLED](#) ()
- [CRoboDacqNet](#) ^ [GetRoboDacq](#) ()
- [CRoboFluidDeviceNet](#) ^ [GetRoboFluidDevice](#) ()
- [CGilsonDeviceNet](#) ^ [GetGilsonDevice](#) ()
- [CMcsBus_ExtensionNet](#) ^ [GetMcsBus_Extension](#) ()

Additional Inherited Members

11.104.1 Detailed Description

[CRoboocyte2DeviceNet](#) is the class to control the MCS Roboocyte2 device

11.104.2 Constructor & Destructor Documentation

11.104.2.1 CRoboocyte2DeviceNet() `CRoboocyte2DeviceNet (void)`

11.104.3 Member Function Documentation

11.104.3.1 GetAxisLED() `bool GetAxisLED ()`

11.104.3.2 GetGilsonDevice() `CGilsonDeviceNet ^ GetGilsonDevice ()`

11.104.3.3 GetMcsBus_Extension() `CMcsBus_ExtensionNet ^ GetMcsBus_Extension ()`


11.104.3.4 GetRoboDacq() `CRoboDacqNet ^ GetRoboDacq ()`

11.104.3.5 GetRoboFluidDevice() `CRoboFluidDeviceNet ^ GetRoboFluidDevice ()`

11.104.3.6 SetAxisLED() `void SetAxisLED (bool onoff)`

11.105 CRoboStatorDeviceNet Class Reference

Inheritance diagram for CRoboStatorDeviceNet:

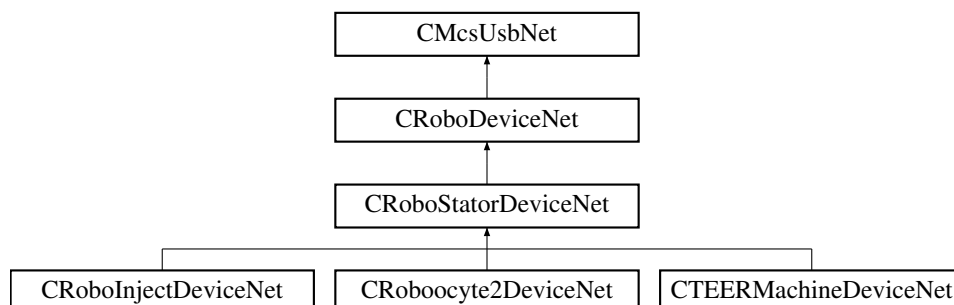

Classes

- class [RoboMainStatorLowLevelCommands](#)

Public Member Functions

- [CRoboStatorDeviceNet](#) (void)
- void [FindReferenceXY](#) ()
- void [FindReferenceXY](#) (int timeout)
- void [FindReferenceZ](#) ()
- void [FindReferenceZ](#) (int timeout)
- void [FindReferenceI](#) ()
- void [FindReferenceI](#) (int timeout)
- unsigned char [HasRefXY](#) ()
- unsigned char [HasRefZ](#) ()
- unsigned char [HasRefI](#) ()
- void [MoveAbsXY](#) (int x, int y)
- void [MoveAbsXY](#) (int x, int y, int timeout)
- void [MoveAbsZ](#) (int z)
- void [MoveAbsZ](#) (int z, int timeout)
- void [MoveAbsI](#) (int i)
- void [MoveAbsI](#) (int i, int timeout)
- void [StopMovementXY](#) ()
- void [StopMovementXY](#) (int timeout)
- void [StopMovementZ](#) ()
- void [StopMovementZ](#) (int timeout)
- void [StopMovementI](#) ()
- void [StopMovementI](#) (int timeout)
- void [SetCurrentAndAirXY](#) (unsigned short onoff)
- void [SetCurrentAndAirXY](#) (unsigned short onoff, int timeout)
- void [GetCurrentPositionXY](#) ([System::Runtime::InteropServices::Out]int% x, [System::Runtime::InteropServices::Out]int% y)
- void [GetCurrentPositionZ](#) ([System::Runtime::InteropServices::Out]int% z)
- void [GetCurrentPositionI](#) ([System::Runtime::InteropServices::Out]int% i)
- void [SetVelocityXY](#) (int v)
- void [SetVelocityZ](#) (int v)
- void [SetVelocityI](#) (int v)
- void [SetSpeedXY](#) (int v)
- void [SetSpeedZ](#) (int v)
- void [SetSpeedI](#) (int v)
- void [SetSpeedNativeXY](#) (int v)
- void [SetSpeedNativeZ](#) (int v)
- void [SetSpeedNativeI](#) (int v)
- void [SetAccelerationXY](#) (int a)
- void [SetAccelerationZ](#) (int a)
- void [SetAccelerationI](#) (int a)
- void [SetAccelerationNativeXY](#) (int a)
- void [SetAccelerationNativeZ](#) (int a)
- void [SetAccelerationNativeI](#) (int a)

Properties

- [RoboMainStatorLowLevelCommands](#)[^] [RoboMainStatorLowLevelCommand](#) [get]

Additional Inherited Members

11.105.1 Constructor & Destructor Documentation

11.105.1.1 CRoboStatorDeviceNet() `CRoboStatorDeviceNet (void)`

11.105.2 Member Function Documentation

11.105.2.1 FindReferenceI() [1/2] `void FindReferenceI ()`

11.105.2.2 FindReferenceI() [2/2] `void FindReferenceI (int timeout)`

11.105.2.3 FindReferenceXY() [1/2] `void FindReferenceXY ()`

11.105.2.4 FindReferenceXY() [2/2] `void FindReferenceXY (int timeout)`

11.105.2.5 FindReferenceZ() [1/2] `void FindReferenceZ ()`

11.105.2.6 FindReferenceZ() [2/2] `void FindReferenceZ (int timeout)`

11.105.2.7 GetCurrentPositionI() `void GetCurrentPositionI ([System::Runtime::InteropServices::Out] int% i)`

11.105.2.8 GetCurrentPositionXY() void GetCurrentPositionXY (
 [System::Runtime::InteropServices::Out] int% x,
 [System::Runtime::InteropServices::Out] int% y)

11.105.2.9 GetCurrentPositionZ() void GetCurrentPositionZ (
 [System::Runtime::InteropServices::Out] int% z)

11.105.2.10 HasRefI() unsigned char HasRefI ()

11.105.2.11 HasRefXY() unsigned char HasRefXY ()

11.105.2.12 HasRefZ() unsigned char HasRefZ ()

11.105.2.13 MoveAbsI() [1/2] void MoveAbsI (
 int i)

11.105.2.14 MoveAbsI() [2/2] void MoveAbsI (
 int i,
 int timeout)

11.105.2.15 MoveAbsXY() [1/2] void MoveAbsXY (
 int x,
 int y)

11.105.2.16 MoveAbsXY() [2/2] void MoveAbsXY (
 int x,
 int y,
 int timeout)

11.105.2.17 MoveAbsZ() [1/2] void MoveAbsZ (
int z)

11.105.2.18 MoveAbsZ() [2/2] void MoveAbsZ (
int z,
int timeout)

11.105.2.19 SetAccelerationI() void SetAccelerationI (
int a)

11.105.2.20 SetAccelerationNativeI() void SetAccelerationNativeI (
int a)

11.105.2.21 SetAccelerationNativeXY() void SetAccelerationNativeXY (
int a)

11.105.2.22 SetAccelerationNativeZ() void SetAccelerationNativeZ (
int a)

11.105.2.23 SetAccelerationXY() void SetAccelerationXY (
int a)

11.105.2.24 SetAccelerationZ() void SetAccelerationZ (
int a)

11.105.2.25 SetCurrentAndAirXY() [1/2] void SetCurrentAndAirXY (
unsigned short onoff)

11.105.2.26 SetCurrentAndAirXY() [2/2] void SetCurrentAndAirXY (
 unsigned short *onoff*,
 int *timeout*)

11.105.2.27 SetSpeedI() void SetSpeedI (
 int *v*)

11.105.2.28 SetSpeedNativeI() void SetSpeedNativeI (
 int *v*)

11.105.2.29 SetSpeedNativeXY() void SetSpeedNativeXY (
 int *v*)

11.105.2.30 SetSpeedNativeZ() void SetSpeedNativeZ (
 int *v*)

11.105.2.31 SetSpeedXY() void SetSpeedXY (
 int *v*)

11.105.2.32 SetSpeedZ() void SetSpeedZ (
 int *v*)

11.105.2.33 SetVelocityI() void SetVelocityI (
 int *v*)

11.105.2.34 SetVelocityXY() void SetVelocityXY (
 int *v*)

11.105.2.35 SetVelocityZ() `void SetVelocityZ (`
 `int v)`

11.105.2.36 StopMovementI() [1/2] `void StopMovementI ()`

11.105.2.37 StopMovementI() [2/2] `void StopMovementI (`
 `int timeout)`

11.105.2.38 StopMovementXY() [1/2] `void StopMovementXY ()`

11.105.2.39 StopMovementXY() [2/2] `void StopMovementXY (`
 `int timeout)`

11.105.2.40 StopMovementZ() [1/2] `void StopMovementZ ()`


11.105.2.41 StopMovementZ() [2/2] `void StopMovementZ (`
 `int timeout)`

11.105.3 Property Documentation

11.105.3.1 RoboMainStatorLowLevelCommand `RoboMainStatorLowLevelCommands^ RoboMainStator←`
`LowLevelCommand [get]`

11.106 CSafeISDeviceNet Class Reference

Inheritance diagram for CSafeISDeviceNet:

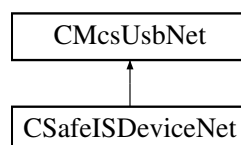

Public Member Functions

- [CSafeISDeviceNet](#) (void)
Initializes a new instance of the [CSafeISDeviceNet](#) class.
- [~CSafeISDeviceNet](#) (void)
Releases unmanaged resources and performs other cleanup operations before the [CSafeISDeviceNet](#) is reclaimed by garbage collection.
- void [SetSwitches](#) (unsigned short switches)
Sets the switches for all electrodes on the device. Do not use during measurement
- void [SetAdcChannels](#) (unsigned char channels)
Sets the ADC channels you want to be sampled
- void [SetAdcSamplePos](#) (array< unsigned short >^ positions)
Sets the sample position of the ADC.
- void [SetDacMode](#) (unsigned char mode)
Sets the DAC mode.
- void [SetDacPulseform](#) (array< short >^ pulseform)
Sets the DAC pulseform.
- void [SetDacPeriode](#) (unsigned int periode)
Sets the DAC periode.

Properties

- [CRoboDeviceNet](#)^ [RoboDevice](#) [get]
Gets the [CRoboDeviceNet](#). Use this to control the syringe.
- [CFluidControlDeviceNet](#)^ [FluidControlDevice](#) [get]
Gets the [CFluidControlDeviceNet](#). Use this to control the valves. Only SetSingleValve is implemented for [CSafeISDeviceNet](#).
- [CMcsUsbDacqNet](#)^ [DacqDevice](#) [get]
Gets the [CMcsUsbDacqNet](#). Use this to control the data aquisition.

Additional Inherited Members

11.106.1 Detailed Description

11.106.2 Constructor & Destructor Documentation

11.106.2.1 [CSafeISDeviceNet](#)() [CSafeISDeviceNet](#) (
void)

Initializes a new instance of the [CSafeISDeviceNet](#) class.

11.106.2.2 [~CSafeISDeviceNet](#)() [~CSafeISDeviceNet](#) (
void)

Releases unmanaged resources and performs other cleanup operations before the [CSafeISDeviceNet](#) is reclaimed by garbage collection.

11.106.3 Member Function Documentation

11.106.3.1 SetAdcChannels() `void SetAdcChannels (`
`unsigned char channels)`

Sets the ADC channels you want to be sampled

Parameters

<i>channels</i>	The bitmap of the 8 channels. Set bit to 1 for the channels you want measure
-----------------	--

11.106.3.2 SetAdcSamplePos() `void SetAdcSamplePos (`
`array< unsigned short >^ positions)`

Sets the sample position of the ADC.

Parameters

<i>positions</i>	The positions in units of 0.1μs.
------------------	----------------------------------

11.106.3.3 SetDacMode() `void SetDacMode (`
`unsigned char mode)`

Sets the DAC mode.

Parameters

<i>mode</i>	The mode: 0 = Impedance ; 1 = Amperometry
-------------	---

11.106.3.4 SetDacPeriode() `void SetDacPeriode (`
`unsigned int periode)`

Sets the DAC periode.

Parameters

<i>periode</i>	The periode in units of 10μs.
----------------	-------------------------------

11.106.3.5 SetDacPulseform() `void SetDacPulseform (array< short >^ pulseform)`

Sets the DAC pulseform.

Parameters

<i>pulseform</i>	The pulseform.
----------------------------------	----------------

11.106.3.6 SetSwitches() `void SetSwitches (unsigned short switches)`

Sets the switches for all electrodes on the device. Do not use during measurement

Parameters

<i>switches</i>	The switches: See Schematics for the meaning
---------------------------------	--

11.106.4 Property Documentation

11.106.4.1 DacqDevice [CMcsUsbDacqNet](#)^ DacqDevice [get]

Gets the [CMcsUsbDacqNet](#). Use this to control the data aquisition.

11.106.4.2 FluidControlDevice [CFluidControlDeviceNet](#)^ FluidControlDevice [get]


Gets the [CFluidControlDeviceNet](#). Use this to control the valves. Only SetSingleValve is implemented for [CSafeISDeviceNet](#).

11.106.4.3 RoboDevice [CRoboDeviceNet](#)^ RoboDevice [get]

Gets the [CRoboDeviceNet](#). Use this to control the syringe.

11.107 CSCUDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CSCUDacqGroupChannelSelectionNet:

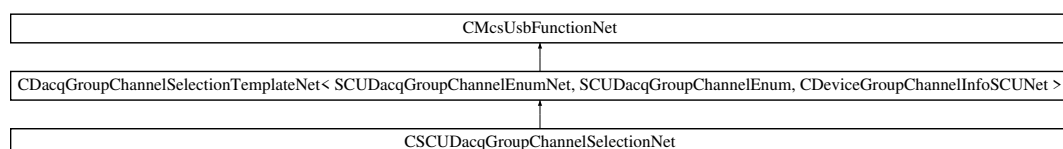

Public Member Functions

- [CSCUDacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)^ mcsusb)

Additional Inherited Members


11.107.1 Constructor & Destructor Documentation

11.107.1.1 [CSCUDacqGroupChannelSelectionNet\(\)](#) [CSCUDacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)^ mcsusb)

11.108 CSCUFunctionNet Class Reference

[CSCUFunctionNet](#) is the class to control the SCU device

Inheritance diagram for [CSCUFunctionNet](#):

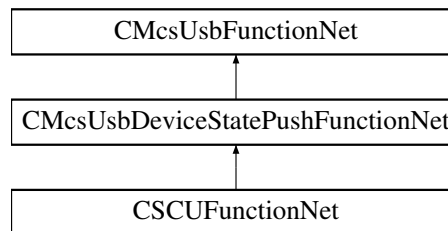

Public Member Functions

- delegate void [OnGetAvailableHeadstages](#) (uint32_t AvailableHeadstages)
- delegate void [OnIsHeadstageAvailable](#) (uint32_t Headstage, bool available)
- [CSCUFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ pSCUFunctionPointerContainer)
 - Initializes a new instance of the [CSCUFunctionNet](#) class.*
- [CSCUFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- virtual [~CSCUFunctionNet](#) ()
- [ICSCUFunctionNet](#) ()
- uint32_t [GetAvailableHeadstages](#) ()
 - Gets a bitmap of available headstages.*
- bool [IsInDacqLegacyMode](#) ()
 - Is the SCU in legacy mode*
- void [SetDacqLegacyMode](#) (bool enable)
 - Enable the SCU legacy mode*
- uint32_t [GetMaxStimulusChannelsPerHeadstage](#) ()
 - Gets the maximal number of stimulation channels a headstage can have.*
- uint32_t [GetMaxNumberOfHeadstages](#) ()
 - Gets the maximal number of headstages.*
- [SCU_HeadstageEnumNet](#) [GetHeadstageID](#) (uint32_t Headstage)

- Gets the headstage fpga ID.*
- bool [IsHeadstageAvailable](#) (uint32_t Headstage)
 - Checks whether the given headstage is available.*
- void [PowerHS](#) (uint32_t Headstage, bool power)
 - Power the HS*
- bool [IsHSPowered](#) (uint32_t Headstage)
 - Is the HS powered*
- bool [HasHSPowerSwitch](#) ()
 - Has SCU HS power switch*
- String ^ [GetHeadstageSerialNumber](#) (uint32_t Headstage)
 - Gets the serial number of a given headstage.*
- uint32_t [GetHeadstageNumberOfAnalogChannels](#) (uint32_t Headstage)
 - Gets the number of analog channels for a given headstage.*
- uint32_t [GetHeadstageNumberOfStimulationChannels](#) (uint32_t Headstage)
 - Gets the number of stimulation channels for a given headstage.*
- uint32_t [GetHeadstageGainInPer mille](#) (uint32_t Headstage)
 - Gets the gain factor in per mille for a given headstage.*
- uint32_t [GetHeadstageAdcRangeInMicroVolt](#) (uint32_t Headstage)
 - Gets the ADC Range in uV for a given headstage.*
- uint32_t [GetHeadstageAdcBits](#) (uint32_t Headstage)
 - Gets the Number of ADC bits for a given headstage.*
- uint32_t [GetHeadstageDacVoltageRangeInMilliVolt](#) (uint32_t Headstage)
 - Gets the DAC Voltage Range in mV for a given headstage.*
- uint32_t [GetHeadstageDacVoltageResolutionInMicroVolt](#) (uint32_t Headstage)
 - Gets the DAC Voltage Resolution in uV for a given headstage.*
- uint32_t [GetHeadstageDacCurrentRangeInMicroAmpere](#) (uint32_t Headstage)
 - Gets the DAC Current Range in uA for a given headstage.*
- uint32_t [GetHeadstageDacCurrentResolutionInNanoAmpere](#) (uint32_t Headstage)
 - Gets the DAC Current Resolution in nA for a given headstage.*
- uint32_t [GetHeadstageDacBits](#) (uint32_t Headstage)
 - Gets the Number of DAC bits for a given headstage.*
- uint32_t [GetHeadstageSamplerate](#) (uint32_t Headstage)
 - Gets the Samplerate of a given headstage.*
- bool [GetHeadstagePowerStateAtStart](#) (uint32_t Headstage)
 - Gets the Power Status at SCU Power on of a given headstage.*
- void [SetHeadstagePowerStateAtStart](#) (uint32_t Headstage, bool Powerstatus)
 - Sets the Power Status at SCU Power on of a given headstage.*
- bool [HasGalvanicIsolation](#) ()
 - Has galvanic isolated hardware*
- bool [HasAnalogOut](#) ()
 - Has AnalogOut hardware*
- void [EnableAnalogOut](#) (bool enable)
 - Enables AnalogOut globally*
- bool [IsAnalogOutEnabled](#) ()
 - Is AnalogOut enabled*
- void [SetAnalogOutDACRange](#) (AnalogOut_DAC_Range_EnumNet range)
 - Sets the analog out DAC range*
- [AnalogOut_DAC_Range_EnumNet](#) [GetAnalogOutDACRange](#) ()
 - Gets the analog out DAC range*
- void [SetAnalogOutADCRange](#) (uint32_t range)
 - Sets the analog out ADC range*

- `uint32_t GetAnalogOutADCRange ()`
Gets the analog out ADC range
- `void AutomaticAnalogOut (bool automatic)`
Sets automatic source channel selection
- `bool IsAutomaticAnalogOut ()`
Is Automatic source channel selection selected
- `void SetAnalogOutChannels (uint32_t out_channel, uint32_t source_channel)`
Set the source channel number for a certain output channel
- `uint32_t GetAnalogOutChannels (uint32_t out_channel)`
Get the connected source channel number for a certain output channel
- `void SetReferenceElectrodeSwitchState (uint32_t Headstage, ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)`
Sets the position of the switch for the reference electrode
- `ReferenceElectrodeSwitchPositionEnumNet GetReferenceElectrodeSwitchState (uint32_t Headstage)`
Gets the position of the switch for the reference electrode
- `void SetReferenceElectrodeMode (uint32_t Headstage, ReferenceElectrodeModeEnumNet NewValue)`
Sets the mode for the reference electrode
- `ReferenceElectrodeModeEnumNet GetReferenceElectrodeMode (uint32_t Headstage)`
Gets the mode for the reference electrode
- `CFilterPropertyNet ^ GetFilterProperty (SCUDacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)`
Gets the filter property
- `array< CFilterPropertyNet ^ > ^ GetFilterProperties (SCUDacqGroupChannelEnumNet GroupID, int filter← Configurations_Length)`
Gets multiple filter properties

Events

- `OnGetAvailableHeadstages ^ GetAvailableHeadstagesEvent [add, remove, raise]`
Event fires when the bitmap of available headstages has changed
- `OnIsHeadstageAvailable ^ IsHeadstageAvailableEvent [add, remove, raise]`
Event fires when 'true' if the headstage is connected for the headstage to query has changed

Additional Inherited Members

11.108.1 Detailed Description

`CSCUFunctionNet` is the class to control the SCU device

11.108.2 Constructor & Destructor Documentation

11.108.2.1 `CSCUFunctionNet()` [1/2] `CSCUFunctionNet (CMcsUsbNet ^ mcsusb, CMcsUsbFunctionPointerContainer ^ pSCUFunctionPointerContainer)`

Initializes a new instance of the `CSCUFunctionNet` class.

11.108.2.2 CSCUFunctionNet() [2/2] `CSCUFunctionNet (
CMcsUsbNet^ mcsusb)`

11.108.2.3 ~CSCUFunctionNet() `virtual ~CSCUFunctionNet () [virtual]`

11.108.2.4 "!CSCUFunctionNet() `!CSCUFunctionNet ()`

11.108.3 Member Function Documentation

11.108.3.1 AutomaticAnalogOut() `void AutomaticAnalogOut (
bool automatic)`

Sets automatic source channel selection

Parameters

<i>automatic</i>	Automatic
------------------	-----------

11.108.3.2 EnableAnalogOut() `void EnableAnalogOut (
bool enable)`

Enables AnalogOut globally

Parameters

<i>enable</i>	Enable
---------------	--------

11.108.3.3 GetAnalogOutADCRange() `uint32_t GetAnalogOutADCRange ()`

Gets the analog out ADC range

Returns

Range

11.108.3.4 GetAnalogOutChannels() `uint32_t GetAnalogOutChannels (`
`uint32_t out_channel)`

Get the connected source channel number for a certain output channel

Parameters

<i>out_channel</i>	Output channel number
--------------------	-----------------------

Returns

Source channel number

11.108.3.5 GetAnalogOutDACRange() `AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ()`

Gets the analog out DAC range

Returns

Range

11.108.3.6 GetAvailableHeadstages() `uint32_t GetAvailableHeadstages ()`

Gets a bitmap of available headstages.

Returns

The bitmap of available headstages.

11.108.3.7 GetFilterProperties() `array<CFilterPropertyNet^> ^ GetFilterProperties (`
`SCUDacqGroupChannelEnumNet GroupID,`
`int filterConfigurations_Length)`

Gets multiple filter properties

Parameters

<i>GroupID</i>	The group ID
<i>filterConfigurations_Length</i>	The maximal length of filterConfigurations.

Returns

array of filter properties

11.108.3.8 GetFilterProperty() `CFilterPropertyNet ^ GetFilterProperty (`
`SCUDacqGroupChannelEnumNet GroupID,`
`uint32_t FilterNumber)`

Gets the filter property

Parameters

<i>GroupID</i>	The group ID
<i>FilterNumber</i>	The filter number

Returns

The filter property

11.108.3.9 GetHeadstageAdcBits() `uint32_t GetHeadstageAdcBits (`
`uint32_t Headstage)`

Gets the Number of ADC bits for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of bits the ADC has for the given headstage.

11.108.3.10 GetHeadstageAdcRangeInMicroVolt() `uint32_t GetHeadstageAdcRangeInMicroVolt (`
`uint32_t Headstage)`

Gets the ADC Range in uV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The ADC Range in uV for the given headstage.

11.108.3.11 GetHeadstageDacBits() `uint32_t GetHeadstageDacBits (`
`uint32_t Headstage)`

Gets the Number of DAC bits for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of bits the DAC has for the given headstage.

11.108.3.12 GetHeadstageDacCurrentRangeInMicroAmpere() `uint32_t GetHeadstageDacCurrentRange↔`
`InMicroAmpere (`
`uint32_t Headstage)`

Gets the DAC Current Range in uA for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Current Range in uA for the given headstage.

11.108.3.13 GetHeadstageDacCurrentResolutionInNanoAmpere() `uint32_t GetHeadstageDacCurrent↔`
`ResolutionInNanoAmpere (`
`uint32_t Headstage)`

Gets the DAC Current Resolution in nA for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Current Resolution in nA for the given headstage.

11.108.3.14 GetHeadstageDacVoltageRangeInMilliVolt() `uint32_t GetHeadstageDacVoltageRangeInMilliVolt (uint32_t Headstage)`

Gets the DAC Voltage Range in mV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Voltage Range in mV for the given headstage.

11.108.3.15 GetHeadstageDacVoltageResolutionInMicroVolt() `uint32_t GetHeadstageDacVoltageResolutionInMicroVolt (uint32_t Headstage)`

Gets the DAC Voltage Resolution in uV for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The DAC Voltage Resolution in uV for the given headstage.

11.108.3.16 GetHeadstageGainInPermille() `uint32_t GetHeadstageGainInPermille (uint32_t Headstage)`

Gets the gain factor in permille for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The gain factor in permille for the given headstage.

11.108.3.17 GetHeadstageID() `SCU_HeadstageIdEnumNet GetHeadstageID (`
`uint32_t Headstage)`

Gets the headstage fpga ID.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The headstage fpga ID.

11.108.3.18 GetHeadstageNumberOfAnalogChannels() `uint32_t GetHeadstageNumberOfAnalogChannels`
`(`
`uint32_t Headstage)`

Gets the number of analog channels for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of analog channels the headstage has.

11.108.3.19 GetHeadstageNumberOfStimulationChannels() `uint32_t GetHeadstageNumberOfStimulation↵`
`Channels (`
`uint32_t Headstage)`

Gets the number of stimulation channels for a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The number of stimulation channels the headstage has.

11.108.3.20 GetHeadstagePowerStateAtStart() `bool GetHeadstagePowerStateAtStart (uint32_t Headstage)`

Gets the Power Status at SCU Power on of a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The Power State at startup for the given headstage: bool false -> off, bool true -> on.

11.108.3.21 GetHeadstageSamplerate() `uint32_t GetHeadstageSamplerate (uint32_t Headstage)`

Gets the Samplerate of a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The samplerate in Hz for the given headstage.

11.108.3.22 GetHeadstageSerialNumber() `String ^ GetHeadstageSerialNumber (uint32_t Headstage)`

Gets the serial number of a given headstage.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

The serial number of the headstage.

11.108.3.23 GetMaxNumberOfHeadstages() `uint32_t GetMaxNumberOfHeadstages ()`

Gets the maximal number of headstages.

Returns

The maximal number of headstages.

11.108.3.24 GetMaxStimulusChannelsPerHeadstage() `uint32_t GetMaxStimulusChannelsPerHeadstage ()`

Gets the maximal number of stimulation channels a headstage can have.

Returns

The maximal number of stimulation channels a headstage can have.

11.108.3.25 GetReferenceElectrodeMode() `ReferenceElectrodeModeEnumNet GetReferenceElectrode↔
Mode (
uint32_t Headstage)`

Gets the mode for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The mode

11.108.3.26 GetReferenceElectrodeSwitchState() `ReferenceElectrodeSwitchPositionEnumNet Get↔
ReferenceElectrodeSwitchState (
uint32_t Headstage)`

Gets the position of the switch for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
------------------	----------------------

Returns

The switch position

11.108.3.27 HasAnalogOut() `bool HasAnalogOut ()`

Has AnalogOut hardware

Returns

Enabled

11.108.3.28 HasGalvanicIsolation() `bool HasGalvanicIsolation ()`

Has galvanic isolated hardware

Returns

Enabled

11.108.3.29 HasHSPowerSwitch() `bool HasHSPowerSwitch ()`

Has SCU HS power switch

Returns

Has Switch

11.108.3.30 IsAnalogOutEnabled() `bool IsAnalogOutEnabled ()`

Is AnalogOut enabled

Returns

Enabled

11.108.3.31 IsAutomaticAnalogOut() `bool IsAutomaticAnalogOut ()`

Is Automatic source channel selection selected

Returns

Automatic

11.108.3.32 IsHeadstageAvailable() `bool IsHeadstageAvailable (
uint32_t Headstage)`

Checks whether the given headstage is available.

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

'true' if the headstage is connected.

11.108.3.33 IsHSPowered() `bool IsHSPowered (`
 `uint32_t Headstage)`

Is the HS powered

Parameters

<i>Headstage</i>	The headstage to query.
------------------	-------------------------

Returns

'true' if the headstage is powered.

11.108.3.34 IsInDacqLegacyMode() `bool IsInDacqLegacyMode ()`

Is the SCU in legacy mode

Returns

Is Enabled

11.108.3.35 OnGetAvailableHeadstages() `delegate void OnGetAvailableHeadstages (`
 `uint32_t AvailableHeadstages)`

11.108.3.36 OnIsHeadstageAvailable() `delegate void OnIsHeadstageAvailable (`
 `uint32_t Headstage,`
 `bool available)`

11.108.3.37 PowerHS() `void PowerHS (`
 `uint32_t Headstage,`
 `bool power)`

Power the HS

Parameters

<i>Headstage</i>	The headstage to query.
<i>power</i>	'true' if the headstage is powered.

11.108.3.38 SetAnalogOutADCRange() `void SetAnalogOutADCRange (`
`uint32_t range)`

Sets the analog out ADC range

Parameters

<i>range</i>	Range
--------------	-------

11.108.3.39 SetAnalogOutChannels() `void SetAnalogOutChannels (`
`uint32_t out_channel,`
`uint32_t source_channel)`

Set the source channel number for a certain output channel

Parameters

<i>out_channel</i>	Output channel number
<i>source_channel</i>	Source channel number

11.108.3.40 SetAnalogOutDACRange() `void SetAnalogOutDACRange (`
`AnalogOut_DAC_Range_EnumNet range)`

Sets the analog out DAC range

Parameters

<i>range</i>	Range
--------------	-------

11.108.3.41 SetDacqLegacyMode() `void SetDacqLegacyMode (`
`bool enable)`

Enable the SCU legacy mode

Parameters

<i>enable</i>	Enable
---------------	--------

11.108.3.42 SetHeadstagePowerStateAtStart() `void SetHeadstagePowerStateAtStart (`
 `uint32_t Headstage,`
 `bool Powerstatus)`

Sets the Power Status at SCU Power on of a given headstage.

Parameters

<i>Headstage</i>	The headstage number
<i>Powerstatus</i>	The Power State at startup for the given headstage: bool false -> off, bool true -> on.

11.108.3.43 SetReferenceElectrodeMode() `void SetReferenceElectrodeMode (`
 `uint32_t Headstage,`
 `ReferenceElectrodeModeEnumNet NewValue)`

Sets the mode for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
<i>NewValue</i>	The mode

11.108.3.44 SetReferenceElectrodeSwitchState() `void SetReferenceElectrodeSwitchState (`
 `uint32_t Headstage,`
 `ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)`

Sets the position of the switch for the reference electrode

Parameters

<i>Headstage</i>	The headstage number
<i>NewSwitchPos</i>	The switch position

11.108.4 Event Documentation

11.108.4.1 GetAvailableHeadstagesEvent `OnGetAvailableHeadstages^ GetAvailableHeadstagesEvent [add], [remove], [raise]`


Event fires when the bitmap of available headstages has changed

11.108.4.2 IsHeadstageAvailableEvent `OnIsHeadstageAvailable^ IsHeadstageAvailableEvent [add], [remove], [raise]`

Event fires when 'true' if the headstage is connected for the headstage to query has changed

11.109 CSerialPortNet Class Reference

Inheritance diagram for CSerialPortNet:

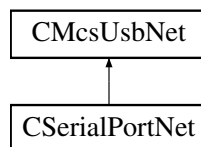

Public Member Functions

- `CSerialPortNet` (void)
- void `Send` (array< byte > ^ buffer)
- void `Send` (String ^ command)
- array< byte > ^ `Receive` (void)
- array< byte > ^ `Receive` (int length)
- String ^ `ReceiveString` (void)
- String ^ `ReceiveString` (int length)
- int `GetBytesAvailable` (void)

Additional Inherited Members

11.109.1 Constructor & Destructor Documentation

11.109.1.1 CSerialPortNet() `CSerialPortNet (void)`

11.109.2 Member Function Documentation

11.109.2.1 GetBytesAvailable() `int GetBytesAvailable (`
`void)`

11.109.2.2 Receive() [1/2] `array<byte> ^ Receive (`
`int length)`

11.109.2.3 Receive() [2/2] `array<byte> ^ Receive (`
`void)`

11.109.2.4 ReceiveString() [1/2] `String ^ ReceiveString (`
`int length)`

11.109.2.5 ReceiveString() [2/2] `String ^ ReceiveString (`
`void)`


11.109.2.6 Send() [1/2] `void Send (`
`array< byte > ^ buffer)`

11.109.2.7 Send() [2/2] `void Send (`
`String ^ command)`

11.110 CStg200xBasicNet Class Reference

Base class for the Stg200x.

Inheritance diagram for CStg200xBasicNet:

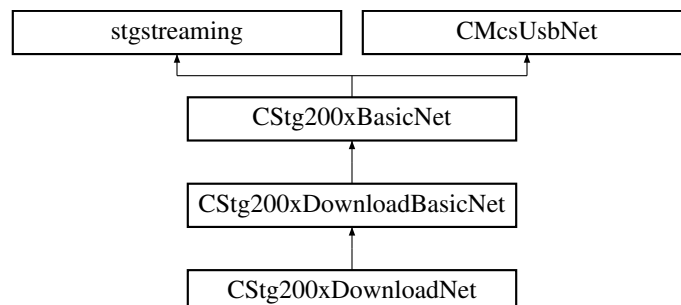

Public Member Functions

- virtual [~CStg200xBasicNet](#) ()
The destructor.
- void [SetOutputRate](#) (uint32_t rate)
Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.
- uint32_t [GetOutputRate](#) ()
Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- void [SendStop](#) (uint32_t triggermap, int options)
Stop some or all triggers of the STG.
- void [GetStgVersionInfo](#) ([Out]String^% SwVersion, [Out]String^% HwVersion)
Queries software and hardware version.
- void [GetAnalogRanges](#) (int channel, [Out]int% URange, [Out]int% IRange)
Gets the range of the analog outputs.
- void [GetAnalogResolution](#) (int channel, [Out]int% URes, [Out]int% IRes)
Gets the resolution of the analog outputs.
- virtual int32_t [GetDACResolution](#) ()
Gets number of bits of the DAC resolution.
- virtual int32_t [GetVoltageRangeInMicroVolt](#) (uint32_t channel)
Gets the Voltage Range of the specified channel in Microvolts.
- virtual int32_t [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)
Gets the Voltage Resolution of the specified channel in Microvolts.
- virtual int32_t [GetCurrentRangeInNanoAmp](#) (uint32_t channel)
Gets the Current Range of the specified channel in Nanoamps.
- virtual int32_t [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)
Gets the Current Resolution of the specified channel in Nanoamps.
- void [GetStgProgramInfo](#) ([Out]bool% IsProgrammed, [Out]System::Runtime::InteropServices::ComTypes::FILETIME% timestamp, [Out]String^% filename, [Out]Guid% guid)
Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.
- void [GetStgProgramInfo](#) ([Out]bool% IsProgrammed, [Out]DateTime% timestamp, [Out]String^% filename, [Out]Guid% guid)
Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.
- void [SetStgProgramInfo](#) (DateTime timestamp, String^ filename, Guid guid)
Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.
- uint32_t [GetAvailableMemory](#) ()
Gets the amount of memory available in the currently selected segment of the STG.
- uint32_t [GetTotalMemory](#) ()
Gets the total amount of memory available on the STG (all segments).
- virtual uint32_t [GetNumberOfAnalogChannels](#) ()
Gets the Number of available analog channels of the device.
- virtual uint32_t [GetNumberOfSyncoutChannels](#) ()
Gets the Number of available syncout channels of the device.
- virtual uint32_t [GetNumberOfTriggerInputs](#) ()
Gets the Number of trigger inputs of the device.
- virtual uint32_t [GetNumberOfHWDACPaths](#) ()

- Gets the Number of HW Stimulation DACs of the device.*

 - virtual uint32_t [GetNumberOfStimulationSourcesPerElectrode](#) ()

Gets the number of stimulation sources (DACs) per electrode.
- virtual void [SetVoltageMode](#) (unsigned int channel)

Sets a channel to voltage mode (STG3008-FA and STG400x only).
- virtual void [SetCurrentMode](#) (unsigned int channel)

Sets a channel to current mode (STG3008-FA and STG400x only).
- virtual void [SetVoltageMode](#) ()

Sets all channels to voltage mode (STG3008-FA and STG400x only).
- virtual void [SetCurrentMode](#) ()

Sets all channels to current mode (STG3008-FA and STG400x only).
- virtual uint32_t [GetCurrentRangeByIndex](#) (uint32_t channel)

Gets the currently selected range index for the current output (STG5 only).
- virtual void [SetCurrentRangeByIndex](#) (uint32_t channel, uint32_t rangeIndex)

Sets the range index for the current output (STG5 only).
- virtual void [SetMeasurementMode](#) (unsigned int channel)

Sets a channel to measurement mode (STG3008-FA).
- virtual void [SetFAAmplification](#) (unsigned int amplification)
- virtual uint32_t [GetFAAmplification](#) ()
- virtual void [SetAutocalibrationDisabled](#) (unsigned int channel, bool disable)

Sets the autocalibration configuration.
- virtual bool [GetAutocalibrationDisabled](#) (unsigned int channel)

Gets the autocalibration configuration.
- virtual void [SetElectrodeMode](#) (uint32_t electrode, array< [ElectrodeModeEnumNet](#) >^ mode)

Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t electrode, [ElectrodeModeEnumNet](#) mode)

Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode, array< [ElectrodeModeEnumNet](#) >^ mode)

Puts an electrode in either automatic or manual mode.
- virtual void [SetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode, [ElectrodeModeEnumNet](#) mode)

Puts an electrode in either automatic or manual mode.
- virtual uint32_t [GetElectrodeMode](#) (uint32_t electrode)

Gets the mode an electrode is in.
- virtual uint32_t [GetElectrodeMode](#) (uint32_t Scu_HS, uint32_t electrode)

Gets the mode an electrode is in.
- virtual void [SetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodelIndex, array< [ElectrodeDacMuxEnumNet](#) >^ dacMux)

Defines the DAC to use for an electrode.
- virtual void [SetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodelIndex, [ElectrodeDacMuxEnumNet](#) dacMux)

Defines the DAC to use for an electrode.
- virtual void [SetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, [ElectrodeDacMuxEnumNet](#) dacMux)

Defines the DAC to use for an electrode.
- virtual void [SetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, array< [ElectrodeDacMuxEnumNet](#) >^ dacMux)

Defines the DAC to use for an electrode.
- virtual [ElectrodeDacMuxEnumNet](#) [GetElectrodeDacMux](#) (uint32_t electrode, uint32_t listmodelIndex)

Gets the DAC which is used for an electrode.
- virtual [ElectrodeDacMuxEnumNet](#) [GetElectrodeDacMux](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex)

- Gets the DAC which is used for an electrode.*
- virtual void [SetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)
Enables or disables the stimulation switch for an electrode.
- virtual void [SetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, bool enable)
Enables or disables the stimulation switch for an electrode.
- virtual void [SetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, bool enable)
Enables or disables the stimulation switch for an electrode.
- virtual void [SetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)
Enables or disables the stimulation switch for an electrode.
- virtual bool [GetElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex)
Gets whether an electrode is enabled or disabled for stimulation.
- virtual bool [GetElectrodeEnable](#) (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodelIndex)
Gets whether an electrode is enabled or disabled for stimulation.
- virtual void [SetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, array< bool >^ enable)
Enables or disables the stimulation switch for an external electrode.
- virtual void [SetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex, bool enable)
Enables or disables the stimulation switch for an external electrode.
- virtual bool [GetExternalElectrodeEnable](#) (uint32_t electrode, uint32_t listmodelIndex)
Gets whether an electrode is enabled or disabled for stimulation.
- virtual void [SetBlankingEnable](#) (uint32_t electrode, bool enable)
Defines whether an electrode should be blanked while stimulation is in progress.
- virtual void [SetBlankingEnable](#) (uint32_t electrode, array< bool >^ enable)
Defines whether an electrode should be blanked while stimulation is in progress.
- virtual void [SetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode, bool enable)
Defines whether an electrode should be blanked while stimulation is in progress.
- virtual void [SetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)
Defines whether an electrode should be blanked while stimulation is in progress.
- virtual bool [GetBlankingEnable](#) (uint32_t electrode)
Gets whether an electrode should be blanked while stimulation is in progress.
- virtual bool [GetBlankingEnable](#) (uint32_t Scu_HS, uint32_t electrode)
Gets whether an electrode should be blanked while stimulation is in progress.
- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t electrode, bool enable)
Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t electrode, array< bool >^ enable)
Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode, bool enable)
Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual void [SetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)
Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual bool [GetEnableAmplifierProtectionSwitch](#) (uint32_t electrode)
Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual bool [GetEnableAmplifierProtectionSwitch](#) (uint32_t Scu_HS, uint32_t electrode)
Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.
- virtual uint32_t [GetNumberOfStimulationElectrodes](#) ()
- template<typename digitalsourceenum >
virtual void [SetTriggerSource](#) (unsigned int triggernum, [DigitalSource](#)< digitalsourceenum >^ triggersource, int bitnum_offset)
- virtual void [SetTriggerSource](#) (unsigned int triggernum, [TriggerSourceEnumNet](#) triggersource, int bitnum_offset)
- virtual void [SetTriggerSource](#) (unsigned int triggernum, [TriggerSourceEnumNet](#) triggersource)

- virtual [TriggerSourceEnumNet GetTriggerSource](#) (unsigned int triggernum)
- virtual void [SetListmodeIndexRange](#) (unsigned int electrodeGroup, unsigned int startIndex, unsigned int endIndex, unsigned int mode)

Define the range of list mode indexes to use for the given electrode group.
- virtual void [GetListmodeIndexRange](#) (unsigned int electrodeGroup, unsigned int &startIndex, unsigned int &endIndex, unsigned int &mode)

Query the range of list mode indexes to use for the given electrode group.
- virtual void [SetListmodeTriggerSource](#) (unsigned int electrodeGroup, [TriggerSourceEnumNet](#) triggersource)

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.
- virtual void [SetListmodeTriggerSource](#) (unsigned int electrodeGroup, [TriggerSourceEnumNet](#) triggersource, int bitnumOffset)

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.
- virtual [TriggerSourceEnumNet GetListmodeTriggerSource](#) (unsigned int electrodeGroup)

Query the currently active signal which triggers the transition from one list mode entry to the next. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.
- virtual void [ListModeSendStart](#) (unsigned int electrodeGroupMask)

Activate (arm) the Listmode for the selected electrode groups.
- virtual void [ListModeSendStop](#) (unsigned int electrodeGroupMask)

Deactivate the Listmode for the selected electrode groups.
- virtual void [SetHeadstage](#) (unsigned int headstage)
- virtual uint32_t [GetHeadstage](#) ()
- virtual void [SetDacAmplificationFactor](#) (uint32_t DacNumber, double Factor)

Set the amplification factor for a DAC.
- virtual double [GetDacAmplificationFactor](#) (uint32_t DacNumber)

Get the amplification factor for a DAC.
- virtual void [SetDigoutMode](#) ([Stg200xDigoutModeEnumNet](#) digoutMode)

Sets the operation mode of the digital output port, can be Monitor, Manual or SyncOut
- virtual [Stg200xDigoutModeEnumNet GetDigoutMode](#) ()

Gets the operation mode of the digital output port, can be Monitor, Manual or SyncOut
- virtual void [SetDigoutValue](#) (uint32_t digoutValue)

Sets the Value on the digital output port when in manual mode.
- virtual uint32_t [GetDigoutValue](#) ()

Gets the Value on the digital output port.
- virtual uint32_t [GetDiginValue](#) ()

Gets the Value on the digital input port.
- virtual void [SetSyncoutMap](#) (uint32_t channel, uint32_t syncoutMap)

Sets the mapping between external syncout outputs and internal syncout channels.
- virtual uint32_t [GetSyncoutMap](#) (uint32_t channel)

Gets the mapping between external syncout outputs and internal syncout channels.

Additional Inherited Members

11.110.1 Detailed Description

Base class for the Stg200x.

From this class all STG related classes are derived: [Mcs.Usb.CStg200xDownloadBasicNet](#) [Mcs.Usb.CStg200xDownloadNet](#) for [Download Mode](#) and [Mcs.Usb.CStg200xStreamingNet](#) for [Streaming Mode](#).

[CStg200xBasicNet](#) is the base class to control MCS STG device.

11.110.2 Constructor & Destructor Documentation

11.110.2.1 `~CStg200xBasicNet()` `virtual ~CStg200xBasicNet () [virtual]`

The destructor.

11.110.3 Member Function Documentation

11.110.3.1 `GetAnalogRanges()` `void GetAnalogRanges (` `int channel,` `[Out] int% URange,` `[Out] int% IRange)`

Gets the range of the analog outputs.

Parameters

<i>channel</i>	The channel which is queried.
<i>URange</i>	The Voltage range in mV.
<i>IRange</i>	The Current range in uA.

11.110.3.2 `GetAnalogResolution()` `void GetAnalogResolution (` `int channel,` `[Out] int% URes,` `[Out] int% IRes)`

Gets the resolution of the analog outputs.

Parameters

<i>channel</i>	The channel which is queried.
<i>URes</i>	The Voltage resolution in mV.
<i>IRes</i>	The Current resolution in uA.

11.110.3.3 `GetAutocalibrationDisabled()` `virtual bool GetAutocalibrationDisabled (` `unsigned int channel) [virtual]`

Gets the autocalibration configuration.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

`true` if autocalibration is disabled.

11.110.3.4 GetAvailableMemory() `uint32_t GetAvailableMemory ()`

Gets the amount of memory available in the currently selected segment of the STG.

Returns

The memory available in the currently selected segment in bytes.

11.110.3.5 GetBlankingEnable() [1/2] `virtual bool GetBlankingEnable (uint32_t electrode) [virtual]`

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

`true` if blanking is enabled while stimulation is in progress.

11.110.3.6 GetBlankingEnable() [2/2] `virtual bool GetBlankingEnable (uint32_t Scu_HS, uint32_t electrode) [virtual]`

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if blanking is enabled while stimulation is in progress.

11.110.3.7 GetCurrentRangeByIndex() `virtual uint32_t GetCurrentRangeByIndex (uint32_t channel) [virtual]`

Gets the currently selected range index for the current output (STG5 only).

Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

Returns

The currently selected range index.

11.110.3.8 GetCurrentRangeInNanoAmp() `virtual int32_t GetCurrentRangeInNanoAmp (uint32_t channel) [virtual]`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.110.3.9 GetCurrentResolutionInNanoAmp() `virtual int32_t GetCurrentResolutionInNanoAmp (uint32_t channel) [virtual]`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.110.3.10 GetDacAmplificationFactor() `virtual double GetDacAmplificationFactor (uint32_t DacNumber) [virtual]`

Get the amplification factor for a DAC.

Parameters

<i>DacNumber</i>	The number of the DAC.
------------------	------------------------

Returns

the amplification factor for the DAC queried, range is from -1.99999 to +1.99999.

11.110.3.11 GetDACResolution() `virtual int32_t GetDACResolution () [virtual]`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.110.3.12 GetDiginValue() `virtual uint32_t GetDiginValue () [virtual]`

Gets the Value on the digital input port.

Returns

The current value on the digital inputs.

11.110.3.13 GetDigoutMode() `virtual Stg200xDigoutModeEnumNet GetDigoutMode () [virtual]`

Gets the operation mode of the digital output port, can be Monitor, Manual or SyncOut

Returns

The current operation mode.

11.110.3.14 GetDigoutValue() `virtual uint32_t GetDigoutValue () [virtual]`

Gets the Value on the digital output port.

Returns

The current value on the digital outputs.

11.110.3.15 GetElectrodeDacMux() [1/2] `virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex) [virtual]`

Gets the DAC which is used for an electrode.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

11.110.3.16 GetElectrodeDacMux() [2/2] `virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex) [virtual]`

Gets the DAC which is used for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

11.110.3.17 GetElectrodeEnable() [1/2] `virtual bool GetElectrodeEnable (`
`uint32_t electrode,`
`uint32_t listmodeIndex) [virtual]`

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

11.110.3.18 GetElectrodeEnable() [2/2] `virtual bool GetElectrodeEnable (`
`uint32_t Scu_HS,`
`uint32_t electrode,`
`uint32_t listmodeIndex) [virtual]`

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

11.110.3.19 GetElectrodeMode() [1/2] `virtual uint32_t GetElectrodeMode (`
`uint32_t electrode) [virtual]`

Gets the mode an electrode is in.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.110.3.20 GetElectrodeMode() [2/2] `virtual uint32_t GetElectrodeMode (`
 `uint32_t Scu_HS,`
 `uint32_t electrode) [virtual]`

Gets the mode an electrode is in.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.110.3.21 GetEnableAmplifierProtectionSwitch() [1/2] `virtual bool GetEnableAmplifierProtection↔`
`Switch (`
 `uint32_t electrode) [virtual]`

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

11.110.3.22 GetEnableAmplifierProtectionSwitch() [2/2] `virtual bool GetEnableAmplifierProtectionSwitch (`
`uint32_t Scu_HS,`
`uint32_t electrode) [virtual]`

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

11.110.3.23 GetExternalElectrodeEnable() `virtual bool GetExternalElectrodeEnable (`
`uint32_t electrode,`
`uint32_t listmodeIndex) [virtual]`

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

<i>electrode</i>	The electrode number.
<i>listmodeIndex</i>	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

11.110.3.24 GetFAAmplification() `virtual uint32_t GetFAAmplification () [virtual]`

11.110.3.25 GetHeadstage() `virtual uint32_t GetHeadstage () [virtual]`

11.110.3.26 GetListModelIndexRange() `virtual void GetListModelIndexRange (unsigned int electrodeGroup, unsigned int & startIndex, unsigned int & endIndex, unsigned int & mode) [virtual]`

Query the range of list mode indexes to use for the given electrode group.

Parameters

<i>electrodeGroup</i>	The electrodegroup for which the range is queried.
-----------------------	--

Parameters

<i>startIndex</i>	The index of the first active element in the listmode list.
<i>endIndex</i>	The index of the last active element in the listmode list.
<i>mode</i>	0 for "start with startIndex", 1 for "start with endIndex".

11.110.3.27 GetListModelTriggerSource() `virtual TriggerSourceEnumNet GetListModelTriggerSource (unsigned int electrodeGroup) [virtual]`

Query the currently active signal which triggers the transition from one list mode entry to the next. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

Parameters

<i>electrodeGroup</i>	The electrodegroup for which the triggersource is queried.
-----------------------	--

Returns

One of the possible sources for the transition.

11.110.3.28 GetNumberOfAnalogChannels() `virtual uint32_t GetNumberOfAnalogChannels () [virtual]`

Gets the Number of available analog channels of the device.

Returns

The number of analog channels.

11.110.3.29 GetNumberOfHWDACPaths() `virtual uint32_t GetNumberOfHWDACPaths () [virtual]`

Gets the Number of HW Stimulation DACs of the device.

Returns

The number of independent HW Stimulation outputs.

11.110.3.30 GetNumberOfStimulationElectrodes() `virtual uint32_t GetNumberOfStimulationElectrodes () [virtual]`

11.110.3.31 GetNumberOfStimulationSourcesPerElectrode() `virtual uint32_t GetNumberOfStimulationSourcesPerElectrode () [virtual]`

Gets the number of stimulation sources (DACs) per electrode.

Returns

The number of stimulation sources (DACs) per electrode.

11.110.3.32 GetNumberOfSyncoutChannels() `virtual uint32_t GetNumberOfSyncoutChannels () [virtual]`

Gets the Number of available syncout channels of the device.

Returns

The number of analog channels.

11.110.3.33 GetNumberOfTriggerInputs() `virtual uint32_t GetNumberOfTriggerInputs () [virtual]`

Gets the Number of trigger inputs of the device.

Returns

The number of trigger inputs.

11.110.3.34 GetOutputRate() `uint32_t GetOutputRate ()`

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Returns

Returns the current output rate in Hz.

11.110.3.35 GetStgProgramInfo() [1/2] `void GetStgProgramInfo ([Out] bool% IsProgrammed, [Out] DateTime% timestamp, [Out] String^% filename, [Out] Guid% guid)`

Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.

Parameters

<i>IsProgrammed</i>	Flag wether download information is valid.
<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.
<i>guid</i>	A GUID.

11.110.3.36 GetStgProgramInfo() [2/2] `void GetStgProgramInfo (`
`[Out] bool% IsProgrammed,`
`[Out] System::Runtime::InteropServices::ComTypes::FILETIME% timestamp,`
`[Out] String^% filename,`
`[Out] Guid% guid)`

Queries Download information from the STG. If download information was stored by the use of [SetStgProgramInfo](#), this function can be used to retrieve it.

Parameters

<i>IsProgrammed</i>	Flag wether download information is valid.
<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.

11.110.3.37 GetStgVersionInfo() `void GetStgVersionInfo (`
`[Out] String^% SwVersion,`
`[Out] String^% HwVersion)`

Queries software and hardware version.

Parameters

<i>SwVersion</i>	The current Software Version of the STG.
<i>HwVersion</i>	The Hardware Revision of the STG.

11.110.3.38 GetSyncoutMap() `virtual uint32_t GetSyncoutMap (`
`uint32_t channel) [virtual]`

Gets the mapping between external syncout outputs and internal syncout channels.

Parameters

<i>channel</i>	The external syncout output channel number (zero based).
----------------	--

Returns

The bitmap of internal syncout channels mapped to channel.

11.110.3.39 GetTotalMemory() `uint32_t GetTotalMemory ()`

Gets the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

11.110.3.40 GetTriggerSource() `virtual TriggerSourceEnumNet GetTriggerSource (unsigned int triggernum) [virtual]`**11.110.3.41 GetVoltageRangeInMicroVolt()** `virtual int32_t GetVoltageRangeInMicroVolt (uint32_t channel) [virtual]`

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

11.110.3.42 GetVoltageResolutionInMicroVolt() `virtual int32_t GetVoltageResolutionInMicroVolt (uint32_t channel) [virtual]`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

11.110.3.43 ListModeSendStart() `virtual void ListModeSendStart (`
`unsigned int electrodeGroupMask) [virtual]`

Activate (arm) the Listmode for the selected electrode groups.

Parameters

<i>electrodeGroupMask</i>	The bitmask of electrode groups for which the listmode is activated.
---------------------------	--

11.110.3.44 ListModeSendStop() `virtual void ListModeSendStop (`
`unsigned int electrodeGroupMask) [virtual]`

Deactivate the Listmode for the selected electrode groups.

Parameters

<i>electrodeGroupMask</i>	The bitmask of electrodegroups for which the listmode is deactivated.
---------------------------	---

11.110.3.45 SendStart() `void SendStart (`
`uint32_t triggermap)`

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

11.110.3.46 SendStop() [1/2] `void SendStop (`
`uint32_t triggermap)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

11.110.3.47 SendStop() [2/2] `void SendStop (`
 `uint32_t triggermap,`
 `int options)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
<i>options</i>	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout associated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active..

11.110.3.48 SetAutocalibrationDisabled() `virtual void SetAutocalibrationDisabled (`
 `unsigned int channel,`
 `bool disable) [virtual]`

Sets the autocalibration configuration.

Parameters

<i>channel</i>	The channel number.
<i>disable</i>	true if autocalibration is to be disabled.

11.110.3.49 SetBlankingEnable() [1/4] `virtual void SetBlankingEnable (`
 `uint32_t electrode,`
 `array< bool >^ enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.50 SetBlankingEnable() [2/4] `virtual void SetBlankingEnable (`
 `uint32_t electrode,`
 `bool enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.51 SetBlankingEnable() [3/4] `virtual void SetBlankingEnable (`
 `uint32_t Scu_HS,`
 `uint32_t electrode,`
 `array< bool >^ enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.52 SetBlankingEnable() [4/4] `virtual void SetBlankingEnable (`
 `uint32_t Scu_HS,`
 `uint32_t electrode,`
 `bool enable) [virtual]`

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.53 SetCurrentMode() [1/2] `virtual void SetCurrentMode () [virtual]`

Sets all channels to current mode (STG3008-FA and STG400x only).

11.110.3.54 SetCurrentMode() [2/2] virtual void SetCurrentMode (
 unsigned int *channel*) [virtual]

Sets a channel to current mode (STG3008-FA and STG400x only).

Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

11.110.3.55 SetCurrentRangeByIndex() virtual void SetCurrentRangeByIndex (
 uint32_t *channel*,
 uint32_t *rangeIndex*) [virtual]

Sets the range index for the current output (STG5 only).

Parameters

<i>channel</i>	The channel to change.
<i>rangeIndex</i>	The new range index.

11.110.3.56 SetDacAmplificationFactor() virtual void SetDacAmplificationFactor (
 uint32_t *DacNumber*,
 double *Factor*) [virtual]

Set the amplification factor for a DAC.

Parameters

<i>DacNumber</i>	The number of the DAC.
<i>Factor</i>	the amplification factor for that DAC, range is from -1.99999 to +1.99999.

11.110.3.57 SetDigoutMode() virtual void SetDigoutMode (
 [Stg200xDigoutModeEnumNet](#) *digoutMode*) [virtual]

Sets the operation mode of the digital output port, can be Monitor, Manual or SyncOut

Parameters

<i>digoutMode</i>	The new operation mode.
-------------------	-------------------------

11.110.3.58 SetDigoutValue() `virtual void SetDigoutValue (uint32_t digoutValue) [virtual]`

Sets the Value on the digital output port when in manual mode.

Parameters

<i>digoutValue</i>	The new value on the digital outputs.
--------------------	---------------------------------------

11.110.3.59 SetElectrodeDacMux() [1/4] `virtual void SetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex, array< ElectrodeDacMuxEnumNet >^ dacMux) [virtual]`

Defines the DAC to use for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

11.110.3.60 SetElectrodeDacMux() [2/4] `virtual void SetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex, ElectrodeDacMuxEnumNet dacMux) [virtual]`

Defines the DAC to use for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

11.110.3.61 SetElectrodeDacMux() [3/4] virtual void SetElectrodeDacMux (
uint32_t *Scu_HS*,
uint32_t *electrode*,
uint32_t *listmodeIndex*,
array< [ElectrodeDacMuxEnumNet](#) >^ *dacMux*) [virtual]

Defines the DAC to use for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

11.110.3.62 SetElectrodeDacMux() [4/4] virtual void SetElectrodeDacMux (
uint32_t *Scu_HS*,
uint32_t *electrode*,
uint32_t *listmodeIndex*,
[ElectrodeDacMuxEnumNet](#) *dacMux*) [virtual]

Defines the DAC to use for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>dacMux</i>	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

11.110.3.63 SetElectrodeEnable() [1/4] `virtual void SetElectrodeEnable (`
 `uint32_t electrode,`
 `uint32_t listmodeIndex,`
 `array< bool >^ enable) [virtual]`

Enables or disables the stimulation switch for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.64 SetElectrodeEnable() [2/4] `virtual void SetElectrodeEnable (`
 `uint32_t electrode,`
 `uint32_t listmodeIndex,`
 `bool enable) [virtual]`

Enables or disables the stimulation switch for an electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.65 SetElectrodeEnable() [3/4] `virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, array< bool >^ enable) [virtual]`

Enables or disables the stimulation switch for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.66 SetElectrodeEnable() [4/4] `virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode,`

```
uint32_t listmodeIndex,  
bool enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.67 SetElectrodeMode() [1/4] `virtual void SetElectrodeMode (`
`uint32_t electrode,`
`array< ElectrodeModeEnumNet >^ mode) [virtual]`

Puts an electrode in either automatic or manual mode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.110.3.68 SetElectrodeMode() [2/4] `virtual void SetElectrodeMode (`
`uint32_t electrode,`
`ElectrodeModeEnumNet mode) [virtual]`

Puts an electrode in either automatic or manual mode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>mode</i>	0 for automatic and 3 for manual mode.
-------------	--

11.110.3.69 SetElectrodeMode() [3/4] virtual void SetElectrodeMode (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 array< [ElectrodeModeEnumNet](#) >^ *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.110.3.70 SetElectrodeMode() [4/4] virtual void SetElectrodeMode (
 uint32_t *Scu_HS*,
 uint32_t *electrode*,
 [ElectrodeModeEnumNet](#) *mode*) [virtual]

Puts an electrode in either automatic or manual mode.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>mode</i>	0 for automatic and 3 for manual mode.
-------------	--

11.110.3.71 SetEnableAmplifierProtectionSwitch() [1/4] `virtual void SetEnableAmplifierProtectionSwitch (`
`uint32_t electrode,`
`array< bool >^ enable) [virtual]`

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.72 SetEnableAmplifierProtectionSwitch() [2/4] `virtual void SetEnableAmplifierProtectionSwitch (`
`uint32_t electrode,`
`bool enable) [virtual]`

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.73 SetEnableAmplifierProtectionSwitch() [3/4] virtual void SetEnableAmplifierProtectionSwitch (

```

    uint32_t Scu_HS,
    uint32_t electrode,
    array< bool >^ enable ) [virtual]

```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.74 SetEnableAmplifierProtectionSwitch() [4/4] virtual void SetEnableAmplifierProtectionSwitch (

```

    uint32_t Scu_HS,
    uint32_t electrode,
    bool enable ) [virtual]

```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

<i>Scu_HS</i>	The SCU headstage number.
---------------	---------------------------

Parameters

<i>electrode</i>	The electrode number.
<i>enable</i>	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.110.3.75 SetExternalElectrodeEnable() [1/2] virtual void SetExternalElectrodeEnable (

```

    uint32_t electrode,
    uint32_t listmodeIndex,
    array< bool >^ enable ) [virtual]

```

Enables or disables the stimulation switch for an external electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.76 SetExternalElectrodeEnable() [2/2] `virtual void SetExternalElectrodeEnable (`
 `uint32_t electrode,`
 `uint32_t listmodeIndex,`
 `bool enable) [virtual]`

Enables or disables the stimulation switch for an external electrode.

Parameters

<i>electrode</i>	The electrode number.
------------------	-----------------------

Parameters

<i>listmodeIndex</i>	The index for listmode.
<i>enable</i>	1 to enable the electrode, 0 to disable.

11.110.3.77 SetFAAmplification() `virtual void SetFAAmplification (`
 `unsigned int amplification) [virtual]`

11.110.3.78 SetHeadstage() `virtual void SetHeadstage (`
 `unsigned int headstage) [virtual]`

11.110.3.79 SetListmodeIndexRange() virtual void SetListmodeIndexRange (
 unsigned int *electrodeGroup*,
 unsigned int *startIndex*,
 unsigned int *endIndex*,
 unsigned int *mode*) [virtual]

Define the range of list mode indexes to use for the given electrode group.

Parameters

<i>electrodeGroup</i>	The electrodegroup for which the range is defined.
-----------------------	--

Parameters

<i>startIndex</i>	The index of the first active element in the listmode list.
<i>endIndex</i>	The index of the last active element in the listmode list.
<i>mode</i>	0 for "start with startIndex", 1 for "start with endIndex".

11.110.3.80 SetListmodeTriggerSource() [1/2] virtual void SetListmodeTriggerSource (
 unsigned int *electrodeGroup*,
 TriggerSourceEnumNet *triggersource*) [virtual]

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use *SCUDigitalSourceEnumNet* and cast to *TriggerSourceEnumNet* if working with an SCU device.

Parameters

<i>electrodeGroup</i>	The electrodegroup for which the triggersource is defined.
-----------------------	--

Parameters

<i>triggersource</i>	One of the possible sources for the transition.
----------------------	---

11.110.3.81 SetListmodeTriggerSource() [2/2] virtual void SetListmodeTriggerSource (
 unsigned int *electrodeGroup*,
 TriggerSourceEnumNet *triggersource*,
 int *bitnumOffset*) [virtual]

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use *SCUDigitalSourceEnumNet* and cast to *TriggerSourceEnumNet* if working with an SCU device.

Parameters

<i>electrodeGroup</i>	The electrodegroup for which the triggersource is defined.
-----------------------	--

Parameters

<i>triggersource</i>	One of the possible sources for the transition.
<i>bitnumOffset</i>	Number to add to the numeric value of the <i>triggersource</i> enum.

11.110.3.82 SetMeasurementMode() `virtual void SetMeasurementMode (unsigned int channel) [virtual]`

Sets a channel to measurement mode (STG3008-FA).

Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

11.110.3.83 SetOutputRate() `void SetOutputRate (uint32_t rate)`

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Parameters

<i>rate</i>	The new output rate in Hz.
-------------	----------------------------

11.110.3.84 SetStgProgramInfo() `void SetStgProgramInfo (DateTime timestamp, String^ filename, Guid guid)`

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

Parameters

<i>timestamp</i>	The timestamp of last download.
<i>filename</i>	The filename of the downloaded waveform.

11.110.3.85 SetSyncoutMap() virtual void SetSyncoutMap (
 uint32_t *channel*,
 uint32_t *syncoutMap*) [virtual]

Sets the mapping between external syncout outputs and internal syncout channels.

Parameters

<i>channel</i>	The external syncout output channel number (zero based).
<i>syncoutMap</i>	A bitmap of internal syncout channels to map to channel.

11.110.3.86 SetTriggerSource() [1/3] virtual void SetTriggerSource (
 unsigned int *triggernum*,
 DigitalSource< digitalsourceenum >^ *triggersource*,
 int *bitnum_offset*) [virtual]

11.110.3.87 SetTriggerSource() [2/3] virtual void SetTriggerSource (
 unsigned int *triggernum*,
 TriggerSourceEnumNet *triggersource*) [virtual]

11.110.3.88 SetTriggerSource() [3/3] virtual void SetTriggerSource (
 unsigned int *triggernum*,
 TriggerSourceEnumNet *triggersource*,
 int *bitnum_offset*) [virtual]

11.110.3.89 SetVoltageMode() [1/2] virtual void SetVoltageMode () [virtual]

Sets all channels to voltage mode (STG3008-FA and STG400x only).

11.110.3.90 SetVoltageMode() [2/2] virtual void SetVoltageMode (
 unsigned int *channel*) [virtual]

Sets a channel to voltage mode (STG3008-FA and STG400x only).


Parameters

<i>channel</i>	The channel to change.
----------------	------------------------

11.111 CStg200xDownloadBasicNet Class Reference

[CStg200xDownloadBasicNet](#) is the base class to control the download mode of the MCS STG device.

Inheritance diagram for CStg200xDownloadBasicNet:

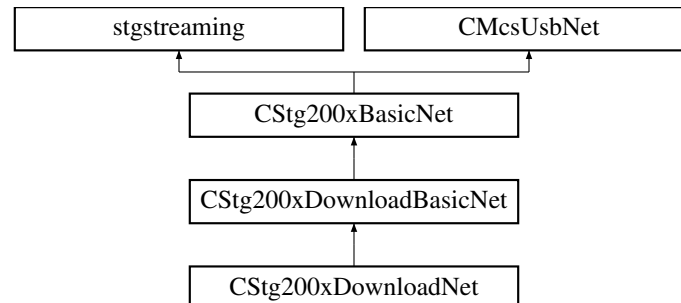

Public Member Functions

- virtual void [SetupTrigger](#) (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)
Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- virtual void [SetupTriggerSingle](#) (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)
Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void [GetTrigger](#) ([Out] array< uint32_t >^% channelmap, [Out] array< uint32_t >^% syncoutmap, [Out] array< uint32_t >^% repeat)
Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void [GetSweepCount](#) ([Out] array< uint32_t >^% sweeps, [Out] array< uint32_t >^% triggers)
Get the sweep and trigger count of the STG.
 - The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
 - The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in [CStg200xDownloadBasicNet::SetupTrigger](#).
- void [ForceStatusEvent](#) ()
Force a status event.
- void [ResetStatus](#) (uint32_t triggermap)
Reset the status flag.
- uint32_t [GetMemoryUsageDAC](#) (uint32_t Channel)
Queries the memory usage of the current segment and selected analog DAC channel.
- uint32_t [GetMemoryUsageSyncout](#) (uint32_t Channel)
Queries the memory usage of the current segment and selected syncout channel.
- virtual void [ClearSyncData](#) (uint32_t channel)
Delete a SyncOut pattern for a channel from STG memory.
- virtual void [SendSyncData](#) (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)
Uploads sync output data to the STG.
Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.
Each datapoint is represented by an integer value and can be either 0 or 1.
The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the [STG200x_SendSyncData32\(\)](#) call instead.
- virtual void [ClearChannelData](#) (uint32_t channel)
Delete a stimulus pattern for a channel from STG memory
- virtual void [SendChannelData](#) (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

- virtual void [EnableAutoReset](#) ()
Enable AutoReset of the STG Status.
- virtual void [DisableAutoReset](#) ()
Disable AutoReset of the STG Status.
- virtual void [SetupRetriggerMode](#) (int8_t trigger, [RetriggerActionEnumNet](#) same_trigger, [RetriggerActionEnumNet](#) other_trigger)
Define the action on triggers while the STG is running.
The STG has three options how to handle a successive trigger while a trigger is active.
 - stop this trigger (default action)
 - restart this trigger
 - ignore the signal
- virtual void [SetupRetriggerMode](#) ([RetriggerActionEnumNet](#) same_trigger, [RetriggerActionEnumNet](#) other_trigger)
Define the action on triggers while the STG is running.
The STG has three options how to handle a successive trigger while a trigger is active.
 - stop this trigger (default action)
 - restart this trigger
 - ignore the signal

Properties

- [CStimulusFunctionNet](#)[^] [Stimulus](#) [get]

Additional Inherited Members

11.111.1 Detailed Description

[CStg200xDownloadBasicNet](#) is the base class to control the download mode of the MCS STG device.

11.111.2 Member Function Documentation

11.111.2.1 ClearChannelData() `virtual void ClearChannelData (uint32_t channel) [virtual]`

Delete a stimulus pattern for a channel from STG memory

Parameters

<i>channel</i>	Specifies the channel to clear.
----------------	---------------------------------

11.111.2.2 ClearSyncData() `virtual void ClearSyncData (
uint32_t channel) [virtual]`

Delete a SyncOut pattern for a channel from STG memory.

Parameters

<i>channel</i>	Specifies the syncout channel to clear.
----------------	---

11.111.2.3 DisableAutoReset() `virtual void DisableAutoReset () [virtual]`

Disable AutoReset of the STG Status.

If autoreset is disabled, the STG status switches to FINISHED after the defined number of sweeps is finished. To switch back to the IDLE status, use CStg200xDownload::ResetStatus()

11.111.2.4 EnableAutoReset() `virtual void EnableAutoReset () [virtual]`

Enable AutoReset of the STG Status.

This is the default on power up. If autoreset is enabled, the STG status switches to FINISHED only for one poll cycle after this, it switches to IDLE automatically.

11.111.2.5 ForceStatusEvent() `void ForceStatusEvent ()`

Force a status event.

Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

11.111.2.6 GetMemoryUsageDAC() `uint32_t GetMemoryUsageDAC (
uint32_t Channel)`

Queries the memory usage of the current segment and selected analog DAC channel.

The currently used memory is reported for the requested channel.

Parameters

<i>Channel</i>	channel for the amount of interested usage.
----------------	---

Returns

Returns the usage in STG memory.

11.111.2.7 GetMemoryUsageSyncout() `uint32_t GetMemoryUsageSyncout (uint32_t Channel)`

Queries the memory usage of the current segment and selected syncout channel.

The currently used memory is reported for the requested channel.

Parameters

<i>Channel</i>	channel for the amount of interested usage.
----------------	---

Returns

Returns the usage in STG memory.

11.111.2.8 GetSweepCount() `void GetSweepCount ([Out] array< uint32_t >^% sweeps, [Out] array< uint32_t >^% triggers)`

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in [CStg200xDownloadBasicNet::SetupTrigger](#).

Parameters

<i>sweeps</i>	on return contains the number of sweeps for each trigger.
<i>triggers</i>	on return contains the number of trigger events seen for each trigger.

11.111.2.9 GetTrigger() `void GetTrigger ([Out] array< uint32_t >^% channelmap, [Out] array< uint32_t >^% syncoutmap, [Out] array< uint32_t >^% repeat)`

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

11.111.2.10 ResetStatus() `void ResetStatus (`
 `uint32_t triggermap)`

Reset the status flag.

Parameters

<i>triggermap</i>	bitmap of trigger for which to reset the status.
-------------------	--

11.111.2.11 SendChannelData() `virtual void SendChannelData (`
 `uint32_t channel,`
 `array< uint16_t >^ pData,`
 `array< uint64_t >^ tData) [virtual]`

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Parameters

<i>channel</i>	Specifies the channel to append the data to.
<i>pData</i>	A list of datapoints.
<i>tData</i>	A list of durations as int64_t. The time is given in units of μ s.

11.111.2.12 SendSyncData() `virtual void SendSyncData (`
 `uint32_t channel,`
 `array< uint16_t >^ pData,`
 `array< uint64_t >^ tData) [virtual]`

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μs . The STG has a resolution of 20 μs . If your application can not handle 64 bit integers, use the STG200x_SendSyncData32() call instead.

Parameters

<i>channel</i>	Specifies the sync output channel to append the data to.
<i>pData</i>	A list of datapoints.
<i>tData</i>	A list of durations as int64_t. The time is given in units of μs .

11.111.2.13 SetupRetriggerMode() [1/2] virtual void SetupRetriggerMode (
 int8_t trigger,
 RetriggerActionEnumNet same_trigger,
 RetriggerActionEnumNet other_trigger) [virtual]

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Parameters

<i>trigger</i>	The trigger to change.
<i>same_trigger</i>	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
<i>other_trigger</i>	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

11.111.2.14 SetupRetriggerMode() [2/2] virtual void SetupRetriggerMode (
 RetriggerActionEnumNet same_trigger,
 RetriggerActionEnumNet other_trigger) [virtual]

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Parameters

<i>same_trigger</i>	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
<i>other_trigger</i>	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

11.111.2.15 SetupTrigger() `virtual void SetupTrigger (`
`uint32_t first_trigger,`
`array< uint32_t >^ channelmap,`
`array< uint32_t >^ syncoutmap,`
`array< uint32_t >^ repeat) [virtual]`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>first_trigger</i>	The number of the first trigger to change.
----------------------	--

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

11.111.2.16 SetupTriggerSingle() `virtual void SetupTriggerSingle (`
`uint32_t trigger,`
`uint32_t channelmap,`
`uint32_t syncoutmap,`
`uint32_t repeat) [virtual]`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>trigger</i>	The trigger to change.
----------------	------------------------

Parameters

<i>channelmap</i>	A bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	A bitmap of syncouts that belong to this trigger.
<i>repeat</i>	The number of times this trigger should be repeated.


11.111.3 Property Documentation

11.111.3.1 Stimulus `CStimulusFunctionNet^` Stimulus [get]

11.112 CStg200xDownloadNet Class Reference

Main class for the STG download mode This class implements the STG download mode interface.

Inheritance diagram for CStg200xDownloadNet:

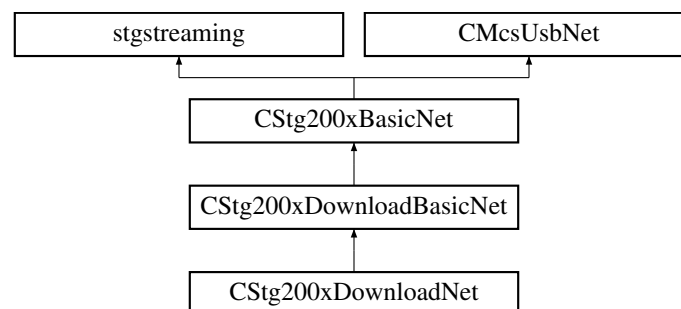

Public Member Functions

- `CStg200xDownloadNet ()`
Use this constructor if you do not want to use the status callback.
- `CStg200xDownloadNet (OnStgPollStatus^ pollStatus)`
Use this constructor if you want to use the status callback.
- `~CStg200xDownloadNet ()`

- void [PrepareAndSendData](#) (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, [STG_DestinationEnumNet](#) destType)
Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.
- void [PrepareAndAppendData](#) (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, [STG_DestinationEnumNet](#) destType)
Prepare and append data to a given channel on the STG.
- void [ClearChannel_PrepareAndSendData](#) (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, [STG_DestinationEnumNet](#) destType, bool doClear)
Prepare and append data to a given channel on the STG.
- void [SendSegmentDefine](#) (array< uint32_t >^ segment_list)
Defines the segment memory layout of the STG.
- void [SendSegmentStart](#) (uint32_t triggermap, uint32_t segment, [Stg200xSegmentFlagsEnumNet](#) segment-flags)
Switchs segment and starts trigger.
- void [SendSegmentSelect](#) (uint32_t segment, [Stg200xSegmentFlagsEnumNet](#) segmentflags)
Switchs segment.
- void [EnableMultiFileMode](#) (uint32_t submode)
Enable the Multi-File mode of the STG.
- void [DisableMultiFileMode](#) ()
Disable the Multi-File mode of the STG
- [StgStatusNet](#) ^ [QueryTriggerstatus](#) ()
- void [SetOutputMap](#) (array< uint32_t >^ ChannelLayout)
- int32_t [GetModuleTemp](#) (unsigned int channel)
- uint32_t [GetModuleCurrent](#) (unsigned int channel)

Events

- [OnStgPollStatus](#)^ [Stg200xPollStatusEvent](#) [add, remove, raise]
- [OnMwPollStatus](#)^ [MwPollStatusEvent](#) [add, remove, raise]

Additional Inherited Members

11.112.1 Detailed Description

Main class for the STG download mode This class implements the STG download mode interface.

11.112.2 Constructor & Destructor Documentation

11.112.2.1 CStg200xDownloadNet() [1/2] [CStg200xDownloadNet](#) ()

Use this constructor if you do not want to use the status callback.

11.112.2.2 CStg200xDownloadNet() [2/2] CStg200xDownloadNet (OnStgPollStatus^ pollStatus)

Use this constructor if you want to use the status callback.

11.112.2.3 ~CStg200xDownloadNet() ~CStg200xDownloadNet ()

11.112.3 Member Function Documentation

11.112.3.1 ClearChannel_PrepareAndSendData() void ClearChannel_PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType, bool doClear)

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μ V, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

11.112.3.2 DisableMultiFileMode() `void DisableMultiFileMode ()`

Disable the Multi-File mode of the STG

Switch the STG back to normal mode. In this mode, trigger inputs are assigned to channels, not to segments.

11.112.3.3 EnableMultiFileMode() `void EnableMultiFileMode (
uint32_t submode)`

Enable the Multi-File mode of the STG.

In Multi-File mode, the trigger inputs switch between segments. To use this mode, define up to as many segments as trigger inputs are available and fill each segment with a stimulus pattern.

Now a trigger on trigger input 1 switches the STG to the first segment and starts all triggers in this segment. Likewise, a trigger on trigger input 2, 3 and 4 selects the respective segment and start all triggers in this segment. So the Multi-File Mode can be used to predefine up to four different stimuli which can be selected without the need for a computer connection.

Parameters

<i>submode</i>	The submode. Submode 0 is regular Multi-File mode as described above, submode 1 is extended Multi-File mode, where the segment is selected based on the digital pattern on the digital inputs. In this mode, 256 different segments can be defined and used.
----------------	--

11.112.3.4 GetModuleCurrent() `uint32_t GetModuleCurrent (
unsigned int channel)`**11.112.3.5 GetModuleTemp()** `int32_t GetModuleTemp (
unsigned int channel)`**11.112.3.6 PrepareAndAppendData()** `void PrepareAndAppendData (
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType)`

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μ V, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

11.112.3.7 PrepareAndSendData() `void PrepareAndSendData (`
 `uint32_t channel,`
 `array< int32_t >^ amplitude,`
 `array< uint64_t >^ duration,`
 `STG_DestinationEnumNet destType)`

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μ V, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μs .
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

11.112.3.8 QueryTriggerstatus() `StgStatusNet ^ QueryTriggerstatus ()`**11.112.3.9 SendSegmentDefine()** `void SendSegmentDefine (array< uint32_t >^ segment_list)`

Defines the segment memory layout of the STG.

On reset, the STG has one segment containing all available memory.

With this command, the STG memory can be devided into several segments. Each segment can be filled with stimulus data.

Parameters

<i>segment_list</i>	The List of memory sizes (one per segment).
---------------------	---

11.112.3.10 SendSegmentSelect() `void SendSegmentSelect (uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)`

Switchs segment.

Parameters

<i>segment</i>	The number of the segment to select.
----------------	--------------------------------------

Parameters

<i>segmentflags</i>	A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.
---------------------	---

11.112.3.11 SendSegmentStart() `void SendSegmentStart (`
 `uint32_t triggermap,`
 `uint32_t segment,`
 `Stg200xSegmentFlagsEnumNet segmentflags)`

Switchs segment and starts trigger.

Parameters

<i>triggermap</i>	A bitmap of triggers that will be started.
-------------------	--

Parameters

<i>segment</i>	The number of the segment to select.
----------------	--------------------------------------

Parameters

<i>segmentflags</i>	A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.
---------------------	---

11.112.3.12 SetOutputMap() `void SetOutputMap (`
 `array< uint32_t >^ ChannelLayout)`


11.112.4 Event Documentation

11.112.4.1 MwPollStatusEvent [OnMwPollStatus](#)[^] MwPollStatusEvent [add], [remove], [raise]

11.112.4.2 Stg200xPollStatusEvent [OnStgPollStatus](#)[^] Stg200xPollStatusEvent [add], [remove], [raise]

11.113 CStimulusFunctionNet Class Reference

Inheritance diagram for CStimulusFunctionNet:

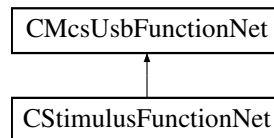

Classes

- class [SidebandData](#)
- class [StimulusDeviceDataAndUnrolledData](#)

Public Member Functions

- [CStimulusFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] stimulusFunctionPointerContainer)
- [CStimulusFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [StartPoll](#) ()
Starts the interrupt fetching thread and delivers events
- void [StopPoll](#) ()
Stops the interrupt fetching thread and delivers events
- void [ForceStatusEvent](#) ()
Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- void [SendStop](#) (uint32_t triggermap, int options)
Stop some or all triggers of the STG.
- void [ClearChannelData](#) (int channel)
Delete a Stimulus Pattern from STG memory
- void [ClearSyncData](#) (int channel)
Delete a Syncout Pattern from STG memory
- void [PrepareAndSendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, [STG_DestinationEnumNet](#) destType)
Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.
- void [PrepareAndAppendData](#) (uint32_t channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, [STG_DestinationEnumNet](#) destType)

Prepare and append data to a given channel on the STG.

- void [ClearChannel_PrepareAndSendData](#) (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, [STG_DestinationEnumNet](#) destType, bool doClear)
- [StimulusDeviceDataAndUnrolledData](#) ^ [PrepareData](#) (int channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, [STG_DestinationEnumNet](#) destType)
- void [SendPreparedData](#) (int channel, [StimulusDeviceDataAndUnrolledData](#)^ device_data_and_unrolled, [STG_DestinationEnumNet](#) destType)
- [SidebandData](#) ^ [CreateSideband](#) (array< int32_t >^ StimulusActive, array< int32_t >^ Syncout, array< uint64_t >^ Duration, uint32_t Bit0Time, uint32_t Bit3Time, uint32_t Bit4Time)

Creates the Sideband Channel for the MEA2100 device.

- void [ClearMultiplexedData](#) ()
- *Clears the Stimulation Memory in the STG device.*
- void [SendMultiplexedData](#) (array< uint16_t >^ data)
- *Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.*
- int [GetMultiplexedDataChannelsInBlock](#) ()
- *Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData*
- int [GetDACResolution](#) ()
- *Gets number of bits of the DAC resolution.*
- int [GetVoltageRangeInMicroVolt](#) (uint32_t channel)
- *Gets the Voltage Range of the specified channel in Microvolts.*
- int [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)
- *Gets the Voltage Resolution of the specified channel in Microvolts.*
- int [GetCurrentRangeInNanoAmp](#) (uint32_t channel)
- *Gets the Current Range of the specified channel in Nanoamps.*
- int [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)
- *Gets the Current Resolution of the specified channel in Nanoamps.*
- void [SetupTrigger](#) (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)
- *Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.*
- void [SetupTriggerSingle](#) (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)
- *Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.*
- uint32_t [GetTotalMemory](#) ()
- *Get the total amount of memory available on the STG (all segments).*
- uint32_t [GetAvailableMemory](#) ()
- *Get the amount of memory available in the currently selected segment of the STG.*
- int [GetNumberOfAnalogChannels](#) ()
- *Get the number of STG channels.*

Events

- [OnStgPollStatus](#)^ [PollStatusEvent](#)

Additional Inherited Members

11.113.1 Constructor & Destructor Documentation

11.113.1.1 CStimulusFunctionNet() [1/2] `CStimulusFunctionNet (`
`CMcsUsbNet^ mcsusb,`
`CMcsUsbFunctionPointerContainer^ stimulusFunctionPointerContainer)`

11.113.1.2 CStimulusFunctionNet() [2/2] `CStimulusFunctionNet (`
`CMcsUsbNet^ mcsusb)`

11.113.2 Member Function Documentation

11.113.2.1 ClearChannel_PrepareAndSendData() `void ClearChannel_PrepareAndSendData (`
`uint32_t channel,`
`array< int32_t >^ amplitude,`
`array< uint64_t >^ duration,`
`STG_DestinationEnumNet destType,`
`bool doClear)`

11.113.2.2 ClearChannelData() `void ClearChannelData (`
`int channel)`

Delete a Stimulus Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

11.113.2.3 ClearMultiplexedData() `void ClearMultiplexedData ()`

Clears the Stimulation Memory in the STG device.

11.113.2.4 ClearSyncData() `void ClearSyncData (`
`int channel)`

Delete a Syncout Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

11.113.2.5 CreateSideband() `SidebandData ^ CreateSideband (`
`array< int32_t >^ StimulusActive,`
`array< int32_t >^ Syncout,`
`array< uint64_t >^ Duration,`
`uint32_t Bit0Time,`
`uint32_t Bit3Time,`
`uint32_t Bit4Time)`

Creates the Sideband Channel for the MEA2100 device.

Each datapoint is represented by an signed 32bit integer value. A value 0 means that the stimulation is active during that time. A value 1 means that the stimulation is not active during that time.

The duration is given as a list of 64 bit integers. Durations are given in units of μs . The STG has a resolution of 20 μs .

Parameters

<i>StimulusActive</i>	A list of datapoints which define weather the Stimulus is active or idle at that time as int32.
-----------------------	---

Parameters

<i>Duration</i>	A list of durations as uint64. The time is given in units of μs .
<i>Bit0Time</i>	Time in μs for which Bit 0 (Blanking) is to be extended.

Parameters

<i>Bit3Time</i>	Time in μs for which Bit 3 (Stimulus Enable) is to be extended.
-----------------	--

Parameters

<i>Bit4Time</i>	Time in μs for which Bit 4 (Stimulus Selector) is to be extended.
-----------------	--

Returns

Error Status. 0 on success.

11.113.2.6 ForceStatusEvent() `void ForceStatusEvent ()`

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

11.113.2.7 GetAvailableMemory() `uint32_t GetAvailableMemory ()`

Get the amount of memory available in the currently selected segment of the STG.

Returns

The total memory available on the STG in bytes.

11.113.2.8 GetCurrentRangeInNanoAmp() `int GetCurrentRangeInNanoAmp (
uint32_t channel)`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.113.2.9 GetCurrentResolutionInNanoAmp() `int GetCurrentResolutionInNanoAmp (
uint32_t channel)`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.113.2.10 GetDACResolution() `int GetDACResolution ()`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.113.2.11 GetMultiplexedDataChannelsInBlock() `int GetMultiplexedDataChannelsInBlock ()`

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in `SendMultiplexedData`

11.113.2.12 GetNumberOfAnalogChannels() `int GetNumberOfAnalogChannels ()`

Get the number of STG channels.

Returns

The number of STG channels.

11.113.2.13 GetTotalMemory() `uint32_t GetTotalMemory ()`

Get the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

11.113.2.14 GetVoltageRangeInMicroVolt() `int GetVoltageRangeInMicroVolt (
uint32_t channel)`

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

11.113.2.15 GetVoltageResolutionInMicroVolt() `int GetVoltageResolutionInMicroVolt (uint32_t channel)`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

11.113.2.16 PrepareAndAppendData() `void PrepareAndAppendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)`

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μ V, thus the possible range is ± 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is ± 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μ V and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μs .
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

11.113.2.17 PrepareAndSendData() `void PrepareAndSendData (`
 `uint32_t channel,`
 `array< int32_t >^ amplitude,`
 `array< uint64_t >^ duration,`
 `STG_DestinationEnumNet destType)`

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 μV , thus the possible range is $\pm 2000 \text{ V}$. When using current stimulation, the values are in multiple of 1 nA, this the possible range is $\pm 2000 \text{ mA}$.

The duration is given as a list of 64 bit integers. Durations are given in units of μs . The STG has a resolution of 20 μs .

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

<i>channel</i>	The channel number to send data to.
----------------	-------------------------------------

Parameters

<i>amplitude</i>	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
------------------	--

Parameters

<i>duration</i>	A list of durations in units of μ s.
<i>destType</i>	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

11.113.2.18 PrepareData() `StimulusDeviceDataAndUnrolledData ^ PrepareData (`
 `int channel,`
 `array< int32_t >^ amplitude,`
 `array< uint64_t >^ duration,`
 `STG_DestinationEnumNet destType)`

11.113.2.19 SendMultiplexedData() `void SendMultiplexedData (`
 `array< uint16_t >^ data)`

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

Parameters

<i>data</i>	Array of data to be sent.
-------------	---------------------------

11.113.2.20 SendPreparedData() `void SendPreparedData (`
 `int channel,`
 `StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled,`
 `STG_DestinationEnumNet destType)`

11.113.2.21 SendStart() `void SendStart (`
 `uint32_t triggermap)`

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

11.113.2.22 SendStop() [1/2] `void SendStop (`
`uint32_t triggermap)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

11.113.2.23 SendStop() [2/2] `void SendStop (`
`uint32_t triggermap,`
`int options)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
<i>options</i>	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout associated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active..

11.113.2.24 SetupTrigger() `void SetupTrigger (`
`uint32_t first_trigger,`
`array< uint32_t >^ channelmap,`
`array< uint32_t >^ syncoutmap,`
`array< uint32_t >^ repeat)`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>first_trigger</i>	The number of the first trigger to change.
----------------------	--

Parameters

<i>channelmap</i>	For each trigger, a bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	For each trigger, a bitmap of syncouts that belong to this trigger.
<i>repeat</i>	For each trigger, define the number of times this trigger should be repeated.

11.113.2.25 SetupTriggerSingle() `void SetupTriggerSingle (`
 `uint32_t trigger,`
 `uint32_t channelmap,`
 `uint32_t syncoutmap,`
 `uint32_t repeat)`

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

<i>trigger</i>	The trigger to change.
----------------	------------------------

Parameters

<i>channelmap</i>	A bitmap of channels that belong to this trigger.
-------------------	---

Parameters

<i>syncoutmap</i>	A bitmap of syncouts that belong to this trigger.
<i>repeat</i>	The number of times this trigger should be repeated.

11.113.2.26 StartPoll() `void StartPoll ()`

Starts the interrupt fetching thread and delivers events


11.113.2.27 StopPoll() `void StopPoll ()`

Stops the interrupt fetching thread and delivers events

11.113.3 Event Documentation**11.113.3.1 PollStatusEvent** `OnStgPollStatus^ PollStatusEvent`**11.114 CSw2to64DeviceNet Class Reference**

The class to control the MCS-USB-Sw2to64 device.

Inheritance diagram for CSw2to64DeviceNet:

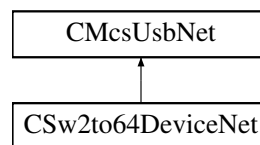**Public Member Functions**

- `CSw2to64DeviceNet ()`
- `~CSw2to64DeviceNet ()`
- unsigned short `GetNumber ()`
Gets the number of channels that can be switched in this box.
- `array< unsigned char > ^ GetChannels ()`
Gets the current switch positions as char array.
- void `SetChannels (array< unsigned char > ^ pattern)`
Sets the switch positions from a char array.
- unsigned char `GetChannel (unsigned short index)`
Gets one current switch position.
- void `SetChannel (unsigned short index, unsigned char pattern)`
Sets one switch position.

Additional Inherited Members**11.114.1 Detailed Description**

The class to control the MCS-USB-Sw2to64 device.

This class controls the settings of the MCS-USB-Sw2to64. The box has two inputs for signals. Each of the 64 outputs can be connected to one of the input signals, could be held open or connected ground. Valid switch states are 0, 1, 2 or 3 for each of the settings.

11.114.2 Constructor & Destructor Documentation

11.114.2.1 CSw2to64DeviceNet() `CSw2to64DeviceNet ()`

11.114.2.2 ~CSw2to64DeviceNet() `~CSw2to64DeviceNet ()`

11.114.3 Member Function Documentation

11.114.3.1 GetChannel() `unsigned char GetChannel (unsigned short index)`

Gets one current switch position.

Parameters

<i>in</i>	<i>index</i>	number of channel to read the switch position from
-----------	--------------	--

Returns

switch position of desired channel

11.114.3.2 GetChannels() `array<unsigned char> ^ GetChannels ()`

Gets the current switch positions as char array.

Returns

array of char with the size of the number of channels, each char has the setting of a channel

11.114.3.3 GetNumber() `unsigned short GetNumber ()`

Gets the number of channels that can be switched in this box.

The box can have a different number of channels it can switch. Up to now usually 64 channels are returned

11.114.3.4 SetChannel() `void SetChannel (unsigned short index, unsigned char pattern)`

Sets one switch position.

Parameters

in	<i>index</i>	number of channel to write the switch position to
in	<i>pattern</i>	switch position of the channel

11.114.3.5 SetChannels() `void SetChannels (`
`array< unsigned char >^ pattern)`

Sets the switch positions from a char array.


Parameters

in	<i>pattern</i>	array of char with the size of the number of channels, each char has the setting of a channel
----	----------------	---

11.115 CTcxDeviceNet Class Reference

Class to control a Temperature Controller (TCX)

Inheritance diagram for CTcxDeviceNet:

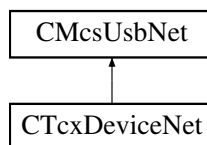

Public Member Functions

- [CTcxDeviceNet \(\)](#)
Initializes a new instance of [CTcxDeviceNet](#) class.
- [~CTcxDeviceNet \(\)](#)
- unsigned int [GetNumControlChannels \(\)](#)
Gets the number of channels the device can control/regulate.
- unsigned int [GetNumMeasureChannels \(\)](#)
Gets the number of channels the device can measure.
- int [GetValue](#) (unsigned int channel)
Gets the temperate of the specified channel in units of 0.1 °C.
- int [GetValueHires](#) (unsigned int channel)
Gets the temperate of the specified channel in units of 0.01 °C.
- int [GetHeaterTemp](#) (unsigned int channel)
Gets the temperate of the specified heater in units of 0.1 °C.
- int [GetHeaterLimit](#) (unsigned int device)
Gets the temperate limit of the specified heater in units of 0.1 °C.
- double [GetMaxHeaterPowerMultiwell \(\)](#)
queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

- void [SetMaxHeaterPowerMultiwell](#) (double MaxPowerWatt)
sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W
- bool [GetHasThermocouple](#) ()
Gets whether the device supports a thermocouple.
- bool [GetEnableHeaterLimit](#) (unsigned int device)
- bool [GetEnableThermocouple](#) (unsigned int device)
- [TcxSensorTypeEnumNet](#) [GetSensorType](#) (unsigned int device)
- String ^ [GetUnit](#) (unsigned int channel)
- unsigned int [GetBoardTemp](#) ()
Gets the temperature of the mainboard in units of 0.1 °C.
- unsigned int [GetVolti](#) (unsigned int channel)
- unsigned int [GetNumDevices](#) ()
- void [SetSetpoint](#) (unsigned int channel, int sp)
Sets the target temperature of specified channel in units of 0.1 °C.
- void [SetDevice](#) (unsigned int channel, int device)
- void [SetOnOff](#) (unsigned int channel, bool on)
Switches the specified channel on or off.
- void [SetCalibration](#) (unsigned int channel, int calib)
- void [SetP](#) (unsigned int device, int p_coeff)
Sets the P-coefficient of the specified device.
- void [SetI](#) (unsigned int device, int i_coeff)
Sets the I-coefficient of the specified device.
- void [SetD](#) (unsigned int device, int d_coeff)
Sets the D-coefficient of the specified device.
- void [SetMaxP](#) (unsigned int device, int maxp)
Sets the maximum heater power of the specified device.
- void [SetHeaterLimit](#) (unsigned int device, int heater_limit)
- void [SetEnableHeaterLimit](#) (unsigned int device, bool enable)
- void [SetEnableThermocouple](#) (unsigned int device, bool enable)
- void [SetSensorType](#) (unsigned int device, [TcxSensorTypeEnumNet](#) type)
- void [SetDevname](#) (unsigned int device, String^ Devicename)
- int [GetSetpoint](#) (unsigned int channel)
Gets the target temperature of specified channel in units of 0.1 °C.
- int [GetDevice](#) (unsigned int channel)
- int [GetOnOff](#) (unsigned int channel)
Gets if the specified channel is on or off.
- int [GetCalibration](#) (unsigned int channel)
- int [GetP](#) (unsigned int device)
Gets the P-coefficient of the specified device.
- int [GetI](#) (unsigned int device)
Gets the I-coefficient of the specified device.
- int [GetD](#) (unsigned int device)
Gets the D-coefficient of the specified device.
- int [GetMaxP](#) (unsigned int device)
Gets the maximum heater power of the specified device.
- String ^ [GetDevname](#) (unsigned int device)
- [TcxDeviceTypeEnumNet](#) [GetDeviceType](#) ()
- int [GetSetpointMin](#) (unsigned int channel)
- int [GetCalibrationMin](#) (unsigned int channel)
- int [GetPMin](#) (unsigned int device)
- int [GetIMin](#) (unsigned int device)
- int [GetDMin](#) (unsigned int device)

- int [GetMaxpMin](#) (unsigned int device)
- int [GetSetpointMax](#) (unsigned int channel)
- int [GetCalibrationMax](#) (unsigned int channel)
- int [GetPMax](#) (unsigned int device)
- int [GetIMax](#) (unsigned int device)
- int [GetDMax](#) (unsigned int device)
- int [GetMaxpMax](#) (unsigned int device)
- int [GetSetpointDecp](#) (unsigned int channel)
- int [GetCalibrationDecp](#) (unsigned int channel)
- int [GetPDecp](#) (unsigned int device)
- int [GetIDecp](#) (unsigned int device)
- int [GetDDecp](#) (unsigned int device)
- int [GetMaxpDecp](#) (unsigned int device)
- int [GetResX](#) (unsigned int channel)
- int [GetResS](#) (unsigned int channel)
- int [GetRes1](#) (unsigned int channel)
- int [GetRes2](#) (unsigned int channel)
- int [GetPwrSet](#) (unsigned int channel)
- int [GetPwrOut](#) (unsigned int channel)
- int [GetDuty](#) (unsigned int channel)
Gets the duty cycle of the heating element.
- int [GetUOut](#) (unsigned int channel)
Gets the voltage on the heating element.
- int [GetIOut](#) (unsigned int channel)
Gets the current through the heating element.
- int [GetROut](#) (unsigned int channel)
Gets the resistance of the heating element.
- int [GetPOut](#) (unsigned int channel)
Gets the output power of the heating element.
- int [GetCurrent](#) (unsigned int channel)
- int [GetThermocoupleTemp](#) (unsigned int channel)
- int [GetThermocoupleTempAbs](#) (unsigned int channel)
- int [GetThermocoupleReferenceTemp](#) (unsigned int channel)
- unsigned int [GetThermocoupleNanovoltPerKelvin](#) (unsigned int channel)
Gets the proportional constant for the thermocouple.
- void [SetThermocoupleNanovoltPerKelvin](#) (unsigned int channel, unsigned int value)
Sets the proportional constant for the thermocouple.
- int [GetThermocoupleCalibration](#) (unsigned int channel)
- void [CalibrateThermocouple](#) (unsigned int channel)
- void [SetDeviceType](#) ([TcxDeviceTypeEnumNet](#) devicetype)
- void [FactoryReset](#) ()

Additional Inherited Members

11.115.1 Detailed Description

Class to control a Temperature Controller (TCX)

11.115.2 Constructor & Destructor Documentation

11.115.2.1 CTcxDeviceNet() `CTcxDeviceNet ()`

Initializes a new instance of `CTcxDeviceNet` class.

11.115.2.2 ~CTcxDeviceNet() `~CTcxDeviceNet ()`**11.115.3 Member Function Documentation****11.115.3.1 CalibrateThermocouple()** `void CalibrateThermocouple (unsigned int channel)`**11.115.3.2 FactoryReset()** `void FactoryReset ()`**11.115.3.3 GetBoardTemp()** `unsigned int GetBoardTemp ()`

Gets the temperate of the mainboard in units of 0.1 °C.

11.115.3.4 GetCalibration() `int GetCalibration (unsigned int channel)`**11.115.3.5 GetCalibrationDecp()** `int GetCalibrationDecp (unsigned int channel)`**11.115.3.6 GetCalibrationMax()** `int GetCalibrationMax (unsigned int channel)`**11.115.3.7 GetCalibrationMin()** `int GetCalibrationMin (unsigned int channel)`

11.115.3.8 GetCurrent() `int GetCurrent (unsigned int channel)`

11.115.3.9 GetD() `int GetD (unsigned int device)`

Gets the D-coefficient of the specified device.

11.115.3.10 GetDDecp() `int GetDDecp (unsigned int device)`

11.115.3.11 GetDevice() `int GetDevice (unsigned int channel)`

11.115.3.12 GetDeviceType() `TcxDeviceTypeEnumNet GetDeviceType ()`

11.115.3.13 GetDevname() `String ^ GetDevname (unsigned int device)`

11.115.3.14 GetDMax() `int GetDMax (unsigned int device)`

11.115.3.15 GetDMin() `int GetDMin (unsigned int device)`

11.115.3.16 GetDuty() `int GetDuty (unsigned int channel)`

Gets the duty cycle of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The duty cycle in percent, the value of $320 * 64$ corresponds to 100 %.

11.115.3.17 GetEnableHeaterLimit() `bool GetEnableHeaterLimit (unsigned int device)`

11.115.3.18 GetEnableThermocouple() `bool GetEnableThermocouple (unsigned int device)`

11.115.3.19 GetHasThermocouple() `bool GetHasThermocouple ()`

Gets whether the device supports a thermocouple.

11.115.3.20 GetHeaterLimit() `int GetHeaterLimit (unsigned int device)`

Gets the temperature limit of the specified heater in units of 0.1 °C.

11.115.3.21 GetHeaterTemp() `int GetHeaterTemp (unsigned int channel)`

Gets the temperature of the specified heater in units of 0.1 °C.

11.115.3.22 GetI() `int GetI (unsigned int device)`

Gets the I-coefficient of the specified device.

11.115.3.23 GetIDecp() `int GetIDecp (`
`unsigned int device)`

11.115.3.24 GetIMax() `int GetIMax (`
`unsigned int device)`

11.115.3.25 GetIMin() `int GetIMin (`
`unsigned int device)`

11.115.3.26 GetIOut() `int GetIOut (`
`unsigned int channel)`

Gets the current through the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The current in units of mA.

11.115.3.27 GetMaxHeaterPowerMultiwell() `double GetMaxHeaterPowerMultiwell ()`

queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

11.115.3.28 GetMaxP() `int GetMaxP (`
`unsigned int device)`

Gets the maximum heater power of the specified device.

11.115.3.29 GetMaxpDecp() `int GetMaxpDecp (`
`unsigned int device)`

11.115.3.30 GetMaxpMax() `int GetMaxpMax (`
`unsigned int device)`

11.115.3.31 GetMaxpMin() `int GetMaxpMin (`
`unsigned int device)`

11.115.3.32 GetNumControlChannels() `unsigned int GetNumControlChannels ()`

Gets the number of channels the device can control/regulate.

11.115.3.33 GetNumDevices() `unsigned int GetNumDevices ()`

11.115.3.34 GetNumMeasureChannels() `unsigned int GetNumMeasureChannels ()`

Gets the number of channels the device can measure.

11.115.3.35 GetOnOff() `int GetOnOff (`
`unsigned int channel)`

Gets if the specified channel is on or off.

11.115.3.36 GetP() `int GetP (`
`unsigned int device)`

Gets the P-coefficient of the specified device.

11.115.3.37 GetPDecp() `int GetPDecp (`
`unsigned int device)`

11.115.3.38 GetPMax() `int GetPMax (`
`unsigned int device)`

11.115.3.39 GetPMin() `int GetPMin (`
`unsigned int device)`

11.115.3.40 GetPOut() `int GetPOut (`
`unsigned int channel)`

Gets the output power of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The resistance in units of mW.

11.115.3.41 GetPwrOut() `int GetPwrOut (unsigned int channel)`

11.115.3.42 GetPwrSet() `int GetPwrSet (unsigned int channel)`

11.115.3.43 GetRes1() `int GetRes1 (unsigned int channel)`

11.115.3.44 GetRes2() `int GetRes2 (unsigned int channel)`

11.115.3.45 GetResS() `int GetResS (unsigned int channel)`

11.115.3.46 GetResX() `int GetResX (unsigned int channel)`

11.115.3.47 GetROut() `int GetROut (unsigned int channel)`

Gets the resistance of the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The resistance in units of 0.1 Ohm.

11.115.3.48 GetSensorType() `TcxSensorTypeEnumNet` GetSensorType (
 unsigned int *device*)

11.115.3.49 GetSetpoint() `int` GetSetpoint (
 unsigned int *channel*)

Gets the target temperate of specified channel in units of 0.1 °C.

11.115.3.50 GetSetpointDecp() `int` GetSetpointDecp (
 unsigned int *channel*)

11.115.3.51 GetSetpointMax() `int` GetSetpointMax (
 unsigned int *channel*)

11.115.3.52 GetSetpointMin() `int` GetSetpointMin (
 unsigned int *channel*)

11.115.3.53 GetThermocoupleCalibration() `int` GetThermocoupleCalibration (
 unsigned int *channel*)

11.115.3.54 GetThermocoupleNanovoltPerKelvin() `unsigned int` GetThermocoupleNanovoltPerKelvin (
 unsigned int *channel*)

Gets the proportional constant for the thermocouple.

Parameters

<i>channel</i>	Thermocouple channel number.
----------------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

11.115.3.55 GetThermocoupleReferenceTemp() `int GetThermocoupleReferenceTemp (unsigned int channel)`

11.115.3.56 GetThermocoupleTemp() `int GetThermocoupleTemp (unsigned int channel)`

11.115.3.57 GetThermocoupleTempAbs() `int GetThermocoupleTempAbs (unsigned int channel)`

11.115.3.58 GetUnit() `String ^ GetUnit (unsigned int channel)`

11.115.3.59 GetUOut() `int GetUOut (unsigned int channel)`

Gets the voltage on the heating element.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

Returns

The voltage in units of mV.

11.115.3.60 GetValue() `int GetValue (`
`unsigned int channel)`

Gets the temperate of the specified channel in units of 0.1 °C.

11.115.3.61 GetValueHires() `int GetValueHires (`
`unsigned int channel)`

Gets the temperate of the specified channel in units of 0.01 °C.

11.115.3.62 GetVolti() `unsigned int GetVolti (`
`unsigned int channel)`

11.115.3.63 SetCalibration() `void SetCalibration (`
`unsigned int channel,`
`int calib)`

11.115.3.64 SetD() `void SetD (`
`unsigned int device,`
`int d_coeff)`

Sets the D-coefficient of the specified device.

11.115.3.65 SetDevice() `void SetDevice (`
`unsigned int channel,`
`int device)`

11.115.3.66 SetDeviceType() `void SetDeviceType (`
`TcxDeviceTypeEnumNet devicetype)`

11.115.3.67 SetDevname() `void SetDevname (`
`unsigned int device,`
`String^ Devicename)`

11.115.3.68 SetEnableHeaterLimit() `void SetEnableHeaterLimit (`
 `unsigned int device,`
 `bool enable)`

11.115.3.69 SetEnableThermocouple() `void SetEnableThermocouple (`
 `unsigned int device,`
 `bool enable)`

11.115.3.70 SetHeaterLimit() `void SetHeaterLimit (`
 `unsigned int device,`
 `int heater_limit)`

11.115.3.71 SetI() `void SetI (`
 `unsigned int device,`
 `int i_coeff)`

Sets the I-coefficient of the specified device.

11.115.3.72 SetMaxHeaterPowerMultiwell() `void SetMaxHeaterPowerMultiwell (`
 `double MaxPowerWatt)`

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

11.115.3.73 SetMaxP() `void SetMaxP (`
 `unsigned int device,`
 `int maxp)`

Sets the maximum heater power of the specified device.

11.115.3.74 SetOnOff() `void SetOnOff (`
 `unsigned int channel,`
 `bool on)`

Switches the specified channel on or off.

Parameters

<i>channel</i>	The channel number.
----------------	---------------------

11.115.3.75 SetP() `void SetP (`
 `unsigned int device,`
 `int p_coeff)`

Sets the P-coefficient of the specified device.

11.115.3.76 SetSensorType() `void SetSensorType (`
 `unsigned int device,`
 `TcxSensorTypeEnumNet type)`

11.115.3.77 SetSetpoint() `void SetSetpoint (`
 `unsigned int channel,`
 `int sp)`

Sets the target temperate of specified channel in units of 0.1 °C.

11.115.3.78 SetThermocoupleNanovoltPerKelvin() `void SetThermocoupleNanovoltPerKelvin (`
 `unsigned int channel,`
 `unsigned int value)`

Sets the proportional constant for the thermocouple.


Parameters

<i>channel</i>	Thermocouple channel number.
<i>value</i>	Proportinal constant in Nanovolt per Kelvin.

11.116 CTEERFunctionNet Class Reference

[CTEERFunctionNet](#) is the class to control the TEER device

Inheritance diagram for CTEERFunctionNet:

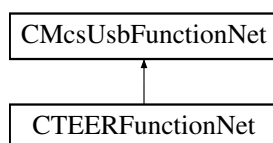

Public Member Functions

- [CTEERFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pTEERFunctionPointer↔
Container)
Initializes a new instance of the [CTEERFunctionNet](#) class.
- [CTEERFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CTEERFunctionNet](#) ()
- [ICTEERFunctionNet](#) ()
- [uint32_t](#) [GetPeriod_us](#) ()
gets the period of TEER stimulation in us
- void [SetPeriod_us](#) ([uint32_t](#) period_us)
sets the period of TEER stimulation in us
- [uint32_t](#) [GetAmplitude_nA](#) ()
gets TEER stimulation amplitude in nA
- void [SetAmplitude_nA](#) ([uint32_t](#) Amplitude_nA)
sets TEER stimulation amplitude in nA
- [TeerWaveformEnumNet](#) [GetWaveform](#) ()
gets TEER stimulation waveform (sine/rect)
- void [SetWaveform](#) ([TeerWaveformEnumNet](#) Waveform)
sets TEER stimulation waveform (sine/rect)
- [TeerClampModeEnumNet](#) [GetClampMode](#) ()
gets TEER clamp mode (voltage/current)
- void [SetClampMode](#) ([TeerClampModeEnumNet](#) ClampMode)
sets TEER clamp mode (voltage/current)
- void [StartSampling](#) ([uint32_t](#) NumberOfCycles)
starts TEER stimulation (duration: n cycles) and samples during last cycle
- void [StopSampling](#) ()
stops TEER stimulation and sampling
- [uint32_t](#) [IsSamplingFinished](#) ()
returns false iff stimulation/sampling is going on, otherwise true
- void [SetControllerParams](#) ([uint32_t](#) P, [uint32_t](#) I, [uint32_t](#) D)
sets PID controller parameters for voltage clamp mode
- void [GetControllerParams](#) ([[System::Runtime::InteropServices::Out](#)][uint32_t](#)% P, [[System::Runtime::InteropServices::Out](#)][uint32_t](#)% I, [[System::Runtime::InteropServices::Out](#)][uint32_t](#)% D)
gets PID controller parameters for voltage clamp mode
- [array< int32_t > ^](#) [GetSampleBufferChunk](#) ([int](#) Buffer_Length)
private function to query max. 100 bytes of sample buffer; called internally
- [array< int32_t > ^](#) [GetSampleVoltageBuffer_uV](#) ([int](#) Buffer_Length)
returns voltage sample buffer (max. 500 values); unit: uV
- [uint32_t](#) [GetMaxChunkSize_Byte](#) ()
private function to be called internally only
- [uint32_t](#) [GetBytesPerSample](#) ()
private function to be called internally only
- [uint32_t](#) [GetNumberOfAvailableSamples](#) ()
private function to be called internally only
- void [SetBufferIndex](#) ([uint32_t](#) NewBufferIndex)
pre-selects sample buffer to be tranferred by [GetSampleVoltageBuffer_uV\(\)](#)
- [uint32_t](#) [GetAdapterCode](#) ()
gets the adapter code
- [uint32_t](#) [GetRotaryPositionCode](#) ()
gets the rotary position code

- void [SetExternalLED](#) (uint32_t NewState)
sets the external LED
- void [SetCurrentEnable](#) (bool NewCurrentEnable)
when disabled, no current will flow through chamber
- bool [GetCurrentEnable](#) ()
when disabled, no current will flow through chamber
- int32_t [GetUptimeSeconds](#) ()
returns time in seconds since device was powered up
- void [StartInternalCalibration](#) ()
starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call
- bool [IsInternalCalibrationFinished](#) ()
queries whether internal calibration has finished
- int [GetDacZero](#) ()
returns DAC-offset (result of internal calibration); use to check for plausibility only
- void [CancelInternalCalibration](#) ()
in case the internal calibration "hangs", this will cancel it
- void [SetLiquidResistance](#) (int32_t NewLiquidResistance_Ohm)
sets the resistance of the liquid in ohms
- int32_t [GetLiquidResistance](#) ()
gets the resistance of the liquid in ohms
- int [GetScaleFactorU1](#) ()
returns U1 scale factor times 10^6 (result of internal calibration)
- int [GetScaleFactorU2](#) ()
returns U2 scale factor times 10^6 (result of internal calibration)
- int [GetAdcOffsetU1](#) ()
returns ADC offset of U1 channel (result of internal calibration)
- int [GetAdcOffsetU2](#) ()
returns ADC offset of U2 channel (result of internal calibration)
- int [GetFrameErrorCounter](#) ()
returns number of times (since bootup) sample memory got overwritten
- int [GetSampleRate](#) ()
returns sample rate in Hz

Additional Inherited Members

11.116.1 Detailed Description

[CTEERFunctionNet](#) is the class to control the TEER device

11.116.2 Constructor & Destructor Documentation

11.116.2.1 CTEERFunctionNet() [1/2] [CTEERFunctionNet](#) (
 CMcsUsbNet^ mcsusb,
 CMcsUsbFunctionPointerContainer^ pTEERFunctionPointerContainer)

Initializes a new instance of the [CTEERFunctionNet](#) class.

11.116.2.2 CTEERFunctionNet() [2/2] `CTEERFunctionNet (CMcsUsbNet^ mcsusb)`

11.116.2.3 ~CTEERFunctionNet() `virtual ~CTEERFunctionNet () [virtual]`

11.116.2.4 !CTEERFunctionNet() `!CTEERFunctionNet ()`

11.116.3 Member Function Documentation

11.116.3.1 CancelInternalCalibration() `void CancelInternalCalibration ()`

in case the internal calibration "hangs", this will cancel it

11.116.3.2 GetAdapterCode() `uint32_t GetAdapterCode ()`

gets the adapter code

Returns

the adapter code

11.116.3.3 GetAdcOffsetU1() `int GetAdcOffsetU1 ()`

returns ADC offset of U1 channel (result of internal calibration)

Returns

the ADC offset for U1

11.116.3.4 GetAdcOffsetU2() `int GetAdcOffsetU2 ()`

returns ADC offset of U2 channel (result of internal calibration)

Returns

the ADC offset for U2

11.116.3.5 GetAmplitude_nA() `uint32_t GetAmplitude_nA ()`

gets TEER stimulation amplitude in nA

Returns

current stimulation amplitude in nA

11.116.3.6 GetBytesPerSample() `uint32_t GetBytesPerSample ()`

private function to be called internally only

Returns

not documented

11.116.3.7 GetClampMode() `TeerClampModeEnumNet GetClampMode ()`

gets TEER clamp mode (voltage/current)

Returns

current TEER clamp mode

11.116.3.8 GetControllerParams() `void GetControllerParams ([System::Runtime::InteropServices::Out] uint32_t% P, [System::Runtime::InteropServices::Out] uint32_t% I, [System::Runtime::InteropServices::Out] uint32_t% D)`

gets PID controller parameters for voltage clamp mode

Parameters

<i>P</i>	the P value
<i>I</i>	the I value
<i>D</i>	the D value

11.116.3.9 GetCurrentEnable() `bool GetCurrentEnable ()`

when disabled, no current will flow through chamber

Returns

false when disabled, true when enabled

11.116.3.10 GetDacZero() `int GetDacZero ()`

returns DAC-offset (result of internal calibration); use to check for plausibility only

Returns

the DAC offset

11.116.3.11 GetFrameErrorCounter() `int GetFrameErrorCounter ()`

returns number of times (since bootup) sample memory got overwritten

Returns

the number of errors

11.116.3.12 GetLiquidResistance() `int32_t GetLiquidResistance ()`

gets the resistance of the liquid in ohms

Returns

the resistance in ohms

11.116.3.13 GetMaxChunkSize_Byte() `uint32_t GetMaxChunkSize_Byte ()`

private function to be called internally only

Returns

not documented

11.116.3.14 GetNumberOfAvailableSamples() `uint32_t GetNumberOfAvailableSamples ()`

private function to be called internally only

Returns

not documented

11.116.3.15 GetPeriod_us() `uint32_t GetPeriod_us ()`

gets the period of TEER stimulation in us

Returns

the period in us

11.116.3.16 GetRotaryPositionCode() `uint32_t GetRotaryPositionCode ()`

gets the rotary position code

Returns

the rotary position code

11.116.3.17 GetSampleBufferChunk() `array<int32_t> ^ GetSampleBufferChunk (
int Buffer_Length)`

private function to query max. 100 bytes of sample buffer; called internally

Parameters

<i>Buffer_Length</i>	The maximal length of Buffer.
----------------------	-------------------------------

Returns

not documented

11.116.3.18 GetSampleRate() `int GetSampleRate ()`

returns sample rate in Hz

Returns

the sample rate in Hz

11.116.3.19 GetSampleVoltageBuffer_uV() `array<int32_t> ^ GetSampleVoltageBuffer_uV (int Buffer_Length)`

returns voltage sample buffer (max. 500 values); unit: uV

Parameters

<i>Buffer_Length</i>	The maximal length of Buffer.
----------------------	-------------------------------

Returns

the voltage sample buffer

11.116.3.20 GetScaleFactorU1() `int GetScaleFactorU1 ()`

returns U1 scale factor times 10^6 (result of internal calibration)

Returns

the U1 scale factor

11.116.3.21 GetScaleFactorU2() `int GetScaleFactorU2 ()`

returns U2 scale factor times 10^6 (result of internal calibration)

Returns

the U2 scale factor

11.116.3.22 GetUptimeSeconds() `int32_t GetUptimeSeconds ()`

returns time in seconds since device was powered up

Returns

seconds since power-on

11.116.3.23 GetWaveform() [TeerWaveformEnumNet](#) GetWaveform ()

gets TEER stimulation waveform (sine/rect)

Returns

waveform enum

11.116.3.24 IsInternalCalibrationFinished() [bool](#) IsInternalCalibrationFinished ()

queries whether internal calibration has finished

Returns

true if calibration has finished

11.116.3.25 IsSamplingFinished() [uint32_t](#) IsSamplingFinished ()

returns false iff stimulation/sampling is going on, otherwise true

Returns

true if sampling is finished

11.116.3.26 SetAmplitude_nA() [void](#) SetAmplitude_nA ([uint32_t](#) *Amplitude_nA*)

sets TEER stimulation amplitude in nA

Parameters

<i>Amplitude_nA</i>	new stimulation amplitude in nA
---------------------	---------------------------------

11.116.3.27 SetBufferIndex() [void](#) SetBufferIndex ([uint32_t](#) *NewBufferIndex*)

pre-selects sample buffer to be transferred by [GetSampleVoltageBuffer_uV\(\)](#)

Parameters

<i>NewBufferIndex</i>	0 - chamber voltage; 1 - compliance voltage
-----------------------	---

11.116.3.28 SetClampMode() `void SetClampMode (`
`TeerClampModeEnumNet ClampMode)`

sets TEER clamp mode (voltage/current)

Parameters

<i>ClampMode</i>	new TEER clamp mode
------------------	---------------------

11.116.3.29 SetControllerParams() `void SetControllerParams (`
`uint32_t P,`
`uint32_t I,`
`uint32_t D)`

sets PID controller parameters for voltage clamp mode

Parameters

<i>P</i>	the P value
<i>I</i>	the I value
<i>D</i>	the D value

11.116.3.30 SetCurrentEnable() `void SetCurrentEnable (`
`bool NewCurrentEnable)`

when disabled, no current will flow through chamber

Parameters

<i>NewCurrentEnable</i>	disabled when false, enabled when true
-------------------------	--

11.116.3.31 SetExternalLED() `void SetExternalLED (`
`uint32_t NewState)`

sets the external LED

Parameters

<i>NewState</i>	state
-----------------	-------

11.116.3.32 SetLiquidResistance() `void SetLiquidResistance (`
`int32_t NewLiquidResistance_Ohm)`

sets the resistance of the liquid in ohms

Parameters

<i>NewLiquidResistance_Ohm</i>	the resistance in ohms
--------------------------------	------------------------

11.116.3.33 SetPeriod_us() `void SetPeriod_us (`
`uint32_t period_us)`

sets the period of TEER stimulation in us

Parameters

<i>period_us</i>	the period in us
------------------	------------------

11.116.3.34 SetWaveform() `void SetWaveform (`
`TeerWaveformEnumNet Waveform)`

sets TEER stimulation waveform (sine/rect)

Parameters

<i>Waveform</i>	waveform enum
-----------------	---------------

11.116.3.35 StartInternalCalibration() `void StartInternalCalibration ()`

starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call

11.116.3.36 StartSampling() `void StartSampling (`
`uint32_t NumberOfCycles)`

starts TEER stimulation (duration: n cycles) and samples during last cycle

Parameters


<i>NumberOfCycles</i>	number of cycles (sine or rect) to output (0 - loop forever)
-----------------------	--

11.116.3.37 StopSampling() void StopSampling ()

stops TEER stimulation and sampling

11.117 CTEERMachineDeviceNet Class Reference

Inheritance diagram for CTEERMachineDeviceNet:

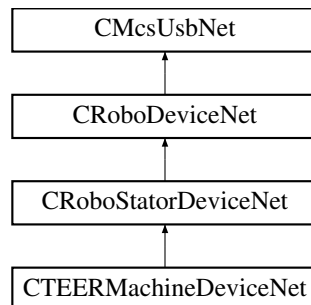**Public Member Functions**

- [CTEERMachineDeviceNet](#) ()
- [~CTEERMachineDeviceNet](#) ()

Properties

- [CTEERFunctionNet](#)[^] [TEERFunctionNet](#) [get]

Additional Inherited Members**11.117.1 Constructor & Destructor Documentation****11.117.1.1 CTEERMachineDeviceNet()** [CTEERMachineDeviceNet](#) ()**11.117.1.2 ~CTEERMachineDeviceNet()** [~CTEERMachineDeviceNet](#) ()


11.117.2 Property Documentation

11.117.2.1 TEERFunctionNet `CTEERFunctionNet^ TEERFunctionNet [get]`

11.118 CUsbDeviceConfigurationFunctionNet Class Reference

`CUsbDeviceConfigurationFunctionNet` is the class to configure the USB firmware

Inheritance diagram for `CUsbDeviceConfigurationFunctionNet`:

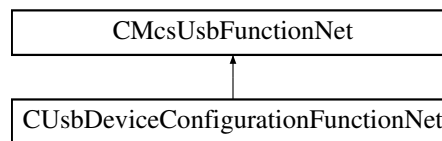

Public Member Functions

- `CUsbDeviceConfigurationFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pUsbDeviceConfigurationFunctionPointerContainer)`
Initializes a new instance of the `CUsbDeviceConfigurationFunctionNet` class.
- `CUsbDeviceConfigurationFunctionNet (CMcsUsbNet^ mcsusb)`
- `virtual ~CUsbDeviceConfigurationFunctionNet ()`
- `!CUsbDeviceConfigurationFunctionNet ()`
- `void SetDeviceName (String^ name)`
sets the USB device name for configurable devices
- `void SetDeviceId (ProductIdEnumNet id)`
sets the USB device name for configurable devices

Additional Inherited Members

11.118.1 Detailed Description

`CUsbDeviceConfigurationFunctionNet` is the class to configure the USB firmware

11.118.2 Constructor & Destructor Documentation

11.118.2.1 CUsbDeviceConfigurationFunctionNet() [1/2] `CUsbDeviceConfigurationFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pUsbDeviceConfigurationFunctionPointerContainer)`

Initializes a new instance of the `CUsbDeviceConfigurationFunctionNet` class.

11.118.2.2 CUsbDeviceConfigurationFunctionNet() [2/2] `CUsbDeviceConfigurationFunctionNet (CMcsUsbNet^ mcsusb)`

11.118.2.3 ~CUsbDeviceConfigurationFunctionNet() `virtual ~CUsbDeviceConfigurationFunctionNet () [virtual]`

11.118.2.4 !CUsbDeviceConfigurationFunctionNet() `!CUsbDeviceConfigurationFunctionNet ()`

11.118.3 Member Function Documentation

11.118.3.1 SetDeviceId() `void SetDeviceId (ProductIdEnumNet id)`

sets the USB device name for configurable devices

Parameters

<i>id</i>	
-----------	--

11.118.3.2 SetDeviceName() `void SetDeviceName (String^ name)`

sets the USB device name for configurable devices


Parameters

<i>name</i>	
-------------	--

11.119 CUsbExceptionNet Class Reference

Exception class that is thrown in case of an USB error.

Inheritance diagram for CUsbExceptionNet:

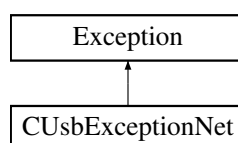

Public Member Functions

- [CUsbExceptionNet](#) (uint32_t status)
Constructor of a CUsbException.
- [CUsbExceptionNet](#) (uint32_t status, String^ message)

Properties

- uint32_t [Status](#) [get]

11.119.1 Detailed Description

Exception class that is thrown in case of an USB error.

11.119.2 Constructor & Destructor Documentation

11.119.2.1 CUsbExceptionNet() [1/2] [CUsbExceptionNet](#) (uint32_t status)

Constructor of a CUsbException.

Parameters

<i>status</i>	the status number
---------------	-------------------

11.119.2.2 CUsbExceptionNet() [2/2] [CUsbExceptionNet](#) (uint32_t status, String^ message)

11.119.3 Property Documentation

11.119.3.1 Status uint32_t Status [get]

11.120 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference

Public Member Functions

- [CVoltageRangeInfoNet](#) (int vr, String^ vrString)

Public Attributes

- int [VoltageRangeInMicroVolt](#)
- String ^ [VoltageRangeDisplayStringMilliVolt](#)

11.120.1 Constructor & Destructor Documentation

11.120.1.1 CVoltageRangeInfoNet() `CVoltageRangeInfoNet (`
 int *vr*,
 String^ *vrString*)


11.120.2 Member Data Documentation

11.120.2.1 VoltageRangeDisplayStringMilliVolt String ^ VoltageRangeDisplayStringMilliVolt

11.120.2.2 VoltageRangeInMicroVolt int VoltageRangeInMicroVolt

11.121 CW2100_FunctionNet Class Reference

Inheritance diagram for CW2100_FunctionNet:

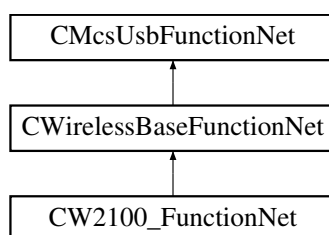**Classes**

- struct [AudioChannelsNet](#)

Public Member Functions

- [CW2100_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ w2100_Function↔
PointerContainer)
- [CW2100_FunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- [array< HeadStageIDType ^ > ^ GetAvailableHeadstages](#) (unsigned int max_length)
- [void SelectHeadstage](#) (unsigned int IDorEntry, int TimeSlotNr)
- [void DeselectHeadstage](#) (int TimeSlotNr)
- [void DeselectAllHeadstages](#) ()
- [HeadStageIDTypeState ^ GetSelectedHeadstageState](#) (int TimeSlotNr)
- [BatteryState ^ GetBatteryState](#) (int TimeSlotNr)
- [System::String ^ GetUserDefinedName](#) (unsigned short ID)
- [System::String ^ GetUserDefinedNameFromSelectedHS](#) (int TimeSlotNr)
- [System::String ^ GetUserDefinedNameCache](#) (unsigned short ID)
- [uint32_t GetUserDefinedNameCache](#) (unsigned short ID, [System::Runtime::InteropServices::Out]System↔
::String^% Name)
- [W2100_StimulusParametersNet ^ GetStiumlusParameters](#) (unsigned short ID)
- [W2100_StimulusParametersNet ^ GetStimulusParametersFromSelectedHS](#) (int TimeSlotNr)
- [W2100_StimulusParametersNet ^ GetStimulusParametersCache](#) (unsigned int typeValue)
- [uint32_t GetStimulusParametersCache](#) (unsigned int typeValue, [System::Runtime::InteropServices::↔
Out]W2100_StimulusParametersNet^% StimulusParameters)
- [void SetSelectedChannels](#) (array< BYTE >^ channels, int TimeSlotNr)
- [array< BYTE > ^ GetSelectedChannels](#) (int TimeSlotNr)
- [void SetMultiHeadstageMode](#) (bool Mode)
- [bool GetMultiHeadstageMode](#) ()
- [void SetHeadstageSamplingActive](#) (bool Active, int TimeSlotNr)
- [bool GetHeadstageSamplingActive](#) (int TimeSlotNr)
- [void SetHeadstageToSleep](#) (unsigned int Sleep16ms, int TimeSlotNr)
- [void SetHeadstageOnOff](#) (unsigned short On, int TimeSlotNr)
- [unsigned short GetHeadstageOnOff](#) (int TimeSlotNr)
- [unsigned int GetAnalogOutChannel](#) ([System::Runtime::InteropServices::Out]int % automatic, unsigned short
index)
- [void SetAnalogOutChannel](#) (int automatic, unsigned short index, unsigned int Channel)
- [array< unsigned int > ^ GetAnalogOutFilter](#) ([System::Runtime::InteropServices::Out]int % automatic)
- [void SetAnalogOutFilter](#) (int automatic, array< unsigned int >^ Coeffs)
- [AnalogOut_DAC_Range_EnumNet GetDacRange](#) ()
- [void SetDacRange](#) (AnalogOut_DAC_Range_EnumNet range)
- [CFilterPropertyNet ^ GetFilterProperty](#) (W2100DacqGroupChannelEnumNet GroupID, unsigned int index)
- [array< CFilterPropertyNet ^ > ^ GetFilterProperties](#) (W2100DacqGroupChannelEnumNet GroupID)
- [void SetAccelGyroEnabled](#) (W2100_Accel_Gyro_Select_EnumNet enable, int TimeSlotNr)
- [W2100_Accel_Gyro_Select_EnumNet GetAccelGyroEnabled](#) (int TimeSlotNr)
- [void SetAccelGyroDesiredRate](#) (int rate, int TimeSlotNr)
- [int GetAccelGyroDesiredRate](#) (int TimeSlotNr)
- [int GetAccelGyroCurrentRate](#) (int TimeSlotNr)
- [void SetAccelRange](#) (int range, int TimeSlotNr)
- [int GetAccelRange](#) (int TimeSlotNr)
- [void SetGyroRange](#) (int range, int TimeSlotNr)
- [int GetGyroRange](#) (int TimeSlotNr)
- [void SetAudioChannels](#) (array< [AudioChannelsNet](#) ^ >^ channels)
- [array< AudioChannelsNet ^ > ^ GetAudioChannels](#) ()
- [unsigned int GetPicFirmwareType](#) (int TimeSlotNr)
- [unsigned int GetFPGA FirmwareType](#) (int TimeSlotNr)

Static Public Member Functions

- static void [ClearUserDefinedNameCache](#) ()
- static void [ClearUserDefinedNameCache](#) (unsigned short ID)
- static void [ClearStimulusParametersCache](#) ()
- static void [ClearStimulusParametersCache](#) (unsigned short ID)

Properties

- [CW2100_StimulatorFunctionNet](#)[^] [Stimulator](#) [get]
- [CPulseGeneratorFunctionNet](#)[^] [PulseGenerator](#) [get]

Additional Inherited Members**11.121.1 Constructor & Destructor Documentation**

11.121.1.1 CW2100_FunctionNet() [1/2] [CW2100_FunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)[^] *w2100_FunctionPointerContainer*)

11.121.1.2 CW2100_FunctionNet() [2/2] [CW2100_FunctionNet](#) (
[CMcsUsbNet](#)[^] *mcsusb*)

11.121.2 Member Function Documentation

11.121.2.1 ClearStimulusParametersCache() [1/2] static void [ClearStimulusParametersCache](#) ()
[static]

11.121.2.2 ClearStimulusParametersCache() [2/2] static void [ClearStimulusParametersCache](#) (
unsigned short *ID*) [static]

11.121.2.3 ClearUserDefinedNameCache() [1/2] static void [ClearUserDefinedNameCache](#) () [static]

11.121.2.4 ClearUserDefinedNameCache() [2/2] static void ClearUserDefinedNameCache (unsigned short *ID*) [static]

11.121.2.5 DeselectAllHeadstages() void DeselectAllHeadstages ()

11.121.2.6 DeselectHeadstage() void DeselectHeadstage (int *TimeSlotNr*)

11.121.2.7 GetAccelGyroCurrentRate() int GetAccelGyroCurrentRate (int *TimeSlotNr*)

11.121.2.8 GetAccelGyroDesiredRate() int GetAccelGyroDesiredRate (int *TimeSlotNr*)

11.121.2.9 GetAccelGyroEnabled() [W2100_Accel_Gyro_Select_EnumNet](#) GetAccelGyroEnabled (int *TimeSlotNr*)

11.121.2.10 GetAccelRange() int GetAccelRange (int *TimeSlotNr*)

11.121.2.11 GetAnalogOutChannel() unsigned int GetAnalogOutChannel ([System::Runtime::InteropServices::Out] int % *automatic*, unsigned short *index*)

11.121.2.12 GetAnalogOutFilter() array<unsigned int> ^ GetAnalogOutFilter ([System::Runtime::InteropServices::Out] int % *automatic*)

11.121.2.13 GetAudioChannels() array<[AudioChannelsNet](#)^> ^ GetAudioChannels ()

11.121.2.14 GetAvailableHeadstages() array<HeadStageIDType^> ^ GetAvailableHeadstages (
 unsigned int *max_length*)

11.121.2.15 GetBatteryState() BatteryState ^ GetBatteryState (
 int *TimeSlotNr*)

11.121.2.16 GetDacRange() AnalogOut_DAC_Range_EnumNet GetDacRange ()

11.121.2.17 GetFilterProperties() array<CFilterPropertyNet^> ^ GetFilterProperties (
 W2100DacqGroupChannelEnumNet *GroupID*)

11.121.2.18 GetFilterProperty() CFilterPropertyNet ^ GetFilterProperty (
 W2100DacqGroupChannelEnumNet *GroupID*,
 unsigned int *index*)

11.121.2.19 GetFPGA FirmwareType() unsigned int GetFPGA FirmwareType (
 int *TimeSlotNr*)

11.121.2.20 GetGyroRange() int GetGyroRange (
 int *TimeSlotNr*)

11.121.2.21 GetHeadstageOnOff() unsigned short GetHeadstageOnOff (
 int *TimeSlotNr*)

11.121.2.22 GetHeadstageSamplingActive() bool GetHeadstageSamplingActive (
 int *TimeSlotNr*)

11.121.2.23 GetMultiHeadstageMode() bool GetMultiHeadstageMode ()

11.121.2.24 GetPicFirmwareType() unsigned int GetPicFirmwareType (
int *TimeSlotNr*)

11.121.2.25 GetSelectedChannels() array<BYTE> ^ GetSelectedChannels (
int *TimeSlotNr*)

11.121.2.26 GetSelectedHeadstageState() HeadStageIDTypeState ^ GetSelectedHeadstageState (
int *TimeSlotNr*)

11.121.2.27 GetStimulusParametersCache() [1/2] W2100_StimulusParametersNet ^ GetStimulus↔
ParametersCache (
unsigned int *typeValue*)

11.121.2.28 GetStimulusParametersCache() [2/2] uint32_t GetStimulusParametersCache (
unsigned int *typeValue*,
[System::Runtime::InteropServices::Out] W2100_StimulusParametersNet^% *Stimulus↔*
Parameters)

11.121.2.29 GetStimulusParametersFromSelectedHS() W2100_StimulusParametersNet ^ GetStimulus↔
ParametersFromSelectedHS (
int *TimeSlotNr*)

11.121.2.30 GetStiumlusParameters() W2100_StimulusParametersNet ^ GetStiumlusParameters (
unsigned short *ID*)

11.121.2.31 GetUserDefinedName() System::String ^ GetUserDefinedName (
unsigned short *ID*)

11.121.2.32 GetUserDefinedNameCache() [1/2] System::String ^ GetUserDefinedNameCache (
unsigned short *ID*)

11.121.2.33 GetUserDefinedNameCache() [2/2] uint32_t GetUserDefinedNameCache (unsigned short *ID*, [System::Runtime::InteropServices::Out] System::String^% *Name*)

11.121.2.34 GetUserDefinedNameFromSelectedHS() System::String ^ GetUserDefinedNameFromSelectedHS (int *TimeSlotNr*)

11.121.2.35 SelectHeadstage() void SelectHeadstage (unsigned int *IDorEntry*, int *TimeSlotNr*)

11.121.2.36 SetAccelGyroDesiredRate() void SetAccelGyroDesiredRate (int *rate*, int *TimeSlotNr*)

11.121.2.37 SetAccelGyroEnabled() void SetAccelGyroEnabled (W2100_Accel_Gyro_Select_EnumNet *enable*, int *TimeSlotNr*)

11.121.2.38 SetAccelRange() void SetAccelRange (int *range*, int *TimeSlotNr*)

11.121.2.39 SetAnalogOutChannel() void SetAnalogOutChannel (int *automatic*, unsigned short *index*, unsigned int *Channel*)

11.121.2.40 SetAnalogOutFilter() void SetAnalogOutFilter (int *automatic*, array< unsigned int >^ *Coeffs*)

11.121.2.41 SetAudioChannels() void SetAudioChannels (
array< AudioChannelsNet[^]>[^] channels)

11.121.2.42 SetDacRange() void SetDacRange (
AnalogOut_DAC_Range_EnumNet range)

11.121.2.43 SetGyroRange() void SetGyroRange (
int range,
int TimeSlotNr)

11.121.2.44 SetHeadstageOnOff() void SetHeadstageOnOff (
unsigned short On,
int TimeSlotNr)

11.121.2.45 SetHeadstageSamplingActive() void SetHeadstageSamplingActive (
bool Active,
int TimeSlotNr)

11.121.2.46 SetHeadstageToSleep() void SetHeadstageToSleep (
unsigned int Sleep16ms,
int TimeSlotNr)

11.121.2.47 SetMultiHeadstageMode() void SetMultiHeadstageMode (
bool Mode)

11.121.2.48 SetSelectedChannels() void SetSelectedChannels (
array< BYTE >[^] channels,
int TimeSlotNr)


11.121.3 Property Documentation

11.121.3.1 PulseGenerator [CPulseGeneratorFunctionNet](#)[^] PulseGenerator [get]

11.121.3.2 Stimulator [CW2100_StimulatorFunctionNet](#)[^] Stimulator [get]

11.122 CW2100_StimulatorFunctionNet Class Reference

Inheritance diagram for CW2100_StimulatorFunctionNet:

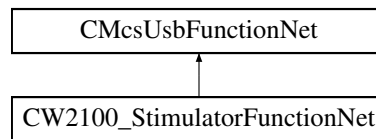

Public Member Functions

- [CW2100_StimulatorFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- void [SendStart](#) (uint32_t triggermap)
Start (Trigger) the STG. The startup delay is in the range of a few ms.
- void [SendStop](#) (uint32_t triggermap)
Stop some or all triggers of the STG.
- [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] [PrepareData](#) (int channel, array< int32_t >[^] amplitude, array< uint64_t >[^] duration, [STG_DestinationEnumNet](#) destType, uint32_t repeat)
- [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] [PrepareDataSync](#) (int channel, array< int32_t >[^] amplitude, array< uint32_t >[^] Sync, array< uint64_t >[^] duration, [STG_DestinationEnumNet](#) destType, uint32_t repeat)
- void [SendPreparedData](#) (int channel, [CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData](#)[^] device_data_and_unrolled, [STG_DestinationEnumNet](#) destType)
- void [ClearChannelData](#) (int channel)
Delete a Stimulus Pattern from STG memory
- int [GetDACResolution](#) ()
Gets number of bits of the DAC resolution.
- int [GetTimeResolutionInNanoSeconds](#) ()
Gets number of bits of the DAC resolution.
- int [GetVoltageRangeInMicroVolt](#) (uint32_t channel)
Gets the Voltage Range of the specified channel in Microvolts.
- int [GetVoltageResolutionInMicroVolt](#) (uint32_t channel)
Gets the Voltage Resolution of the specified channel in Microvolts.
- int [GetCurrentRangeInNanoAmp](#) (uint32_t channel)
Gets the Current Range of the specified channel in Nanoamps.
- int [GetCurrentResolutionInNanoAmp](#) (uint32_t channel)
Gets the Current Resolution of the specified channel in Nanoamps.
- uint32_t [GetNumberOfAnalogChannels](#) ()
- uint32_t [GetNumberOfSyncoutChannels](#) ()
- uint32_t [GetNumberOfTriggerInputs](#) ()
- void [SelectTimeSlot](#) (int TimeSlotNr)
- int [GetTimeSlot](#) ()
- uint32_t [GetStimulationPatternMemory](#) ()

- `uint32_t GetBoostPreTime ()`
- `uint32_t GetBoostAlwaysOnMode ()`
- `void SetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, W2100DigitalSourceEnumNet digstream_source, int bitnumber_offset)`
- `void GetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% digstream_source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)`
- `void SetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, DigitalStimulatorTriggerSlopeEnumNet slope)`
- `DigitalStimulatorTriggerSlopeEnumNet GetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number)`
- `void StartPoll ()`
- `void StopPoll ()`

Static Public Attributes

- `static const uint32_t BOOST_BIT = (1 << 0)`
- `static const uint32_t GND_SWITCH_BIT = (1 << 1)`
- `static const uint32_t SYNC_BIT0 = (1 << 2)`
- `static const uint32_t SYNC_BIT1 = (1 << 3)`

Events

- `OnStgPollStatus^ PollStatusEvent`

Additional Inherited Members

11.122.1 Constructor & Destructor Documentation

11.122.1.1 CW2100_StimulatorFunctionNet() `CW2100_StimulatorFunctionNet (CMcsUsbNet^ mcsusb)`

11.122.2 Member Function Documentation

11.122.2.1 ClearChannelData() `void ClearChannelData (int channel)`

Delete a Stimulus Pattern from STG memory

Parameters

<i>channel</i>	specifies the channel to clear.
----------------	---------------------------------

11.122.2.2 GetBoostAlwaysOnMode() `uint32_t GetBoostAlwaysOnMode ()`

11.122.2.3 GetBoostPreTime() `uint32_t GetBoostPreTime ()`

11.122.2.4 GetCurrentRangeInNanoAmp() `int GetCurrentRangeInNanoAmp (
uint32_t channel)`

Gets the Current Range of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.122.2.5 GetCurrentResolutionInNanoAmp() `int GetCurrentResolutionInNanoAmp (
uint32_t channel)`

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.122.2.6 GetDACResolution() `int GetDACResolution ()`

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.122.2.7 GetDigitalStimulatorTrigger() void GetDigitalStimulatorTrigger (
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
[System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% digstream_↵
source,
[System::Runtime::InteropServices::Out] int% bitnumber_offset)

11.122.2.8 GetDigitalStimulatorTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetDigital↵
StimulatorTriggerSlope (
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number)

11.122.2.9 GetNumberOfAnalogChannels() uint32_t GetNumberOfAnalogChannels ()

11.122.2.10 GetNumberOfSyncoutChannels() uint32_t GetNumberOfSyncoutChannels ()

11.122.2.11 GetNumberOfTriggerInputs() uint32_t GetNumberOfTriggerInputs ()

11.122.2.12 GetStimulationPatternMemory() uint32_t GetStimulationPatternMemory ()

11.122.2.13 GetTimeResolutionInNanoSeconds() int GetTimeResolutionInNanoSeconds ()

Gets number of bits of the DAC resolution.

Returns

The time resolution in ns.

11.122.2.14 GetTimeSlot() int GetTimeSlot ()

11.122.2.15 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (
uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

11.122.2.16 GetVoltageResolutionInMicroVolt() `int GetVoltageResolutionInMicroVolt (uint32_t channel)`

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

<i>channel</i>	Channel which is queried.
----------------	---------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

11.122.2.17 PrepareData() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareData (`

```
int channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType,
uint32_t repeat )
```

11.122.2.18 PrepareDataSync() `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareDataSync (`

```
int channel,
array< int32_t >^ amplitude,
array< uint32_t >^ Sync,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType,
uint32_t repeat )
```

11.122.2.19 SelectTimeSlot() `void SelectTimeSlot (int TimeSlotNr)`

11.122.2.20 SendPreparedData() `void SendPreparedData (`
 `int channel,`
 `CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled,`
 `STG_DestinationEnumNet destType)`

11.122.2.21 SendStart() `void SendStart (`
 `uint32_t triggermap)`

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be started.
-------------------	---

11.122.2.22 SendStop() `void SendStop (`
 `uint32_t triggermap)`

Stop some or all triggers of the STG.

Parameters

<i>triggermap</i>	A bitmap of triggers which will be stopped.
-------------------	---

11.122.2.23 SetDigitalStimulatorTrigger() `void SetDigitalStimulatorTrigger (`
 `int TimeSlotNr,`
 `DigitalStimulatorTriggerEventEnumNet trigger_event,`
 `int trigger_number,`
 `W2100DigitalSourceEnumNet digstream_source,`
 `int bitnumber_offset)`

11.122.2.24 SetDigitalStimulatorTriggerSlope() `void SetDigitalStimulatorTriggerSlope (`
 `int TimeSlotNr,`
 `DigitalStimulatorTriggerEventEnumNet trigger_event,`
 `int trigger_number,`
 `DigitalStimulatorTriggerSlopeEnumNet slope)`

11.122.2.25 StartPoll() `void StartPoll ()`

11.122.2.26 StopPoll() `void StopPoll ()`

11.122.3 Member Data Documentation

11.122.3.1 BOOST_BIT `const uint32_t BOOST_BIT = (1 << 0) [static]`

11.122.3.2 GND_SWITCH_BIT `const uint32_t GND_SWITCH_BIT = (1 << 1) [static]`

11.122.3.3 SYNC_BIT0 `const uint32_t SYNC_BIT0 = (1 << 2) [static]`


11.122.3.4 SYNC_BIT1 `const uint32_t SYNC_BIT1 = (1 << 3) [static]`

11.122.4 Event Documentation

11.122.4.1 PollStatusEvent `OnStgPollStatus^ PollStatusEvent`

11.123 CW2100DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CW2100DacqGroupChannelSelectionNet:

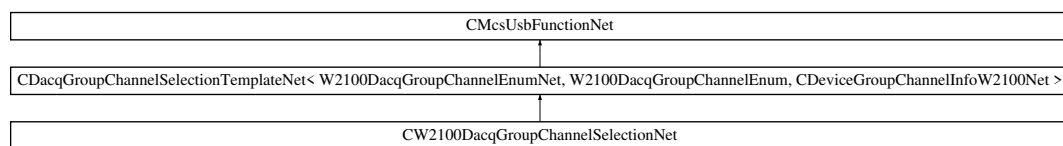

Public Member Functions

- [CW2100DacqGroupChannelSelectionNet](#) ([CMcsUsbNet](#)[^] mcsusb)

Additional Inherited Members


11.123.1 Constructor & Destructor Documentation

11.123.1.1 CW2100DacqGroupChannelSelectionNet() `CW2100DacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)`

11.124 CWarnerUssingDeviceNet Class Reference

`CWarnerUssingDeviceNet` is the class to control the Ussing device

Inheritance diagram for `CWarnerUssingDeviceNet`:

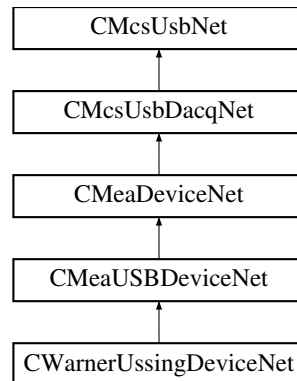

Public Member Functions

- `CWarnerUssingDeviceNet ()`
Initializes a new instance of the `CWarnerUssingDeviceNet` class.
- virtual `~CWarnerUssingDeviceNet ()`
- `ICWarnerUssingDeviceNet ()`

Properties

- `CWarnerUssingFunctionNet^ WarnerUssingFunction` [get]

Additional Inherited Members

11.124.1 Detailed Description

`CWarnerUssingDeviceNet` is the class to control the Ussing device

11.124.2 Constructor & Destructor Documentation

11.124.2.1 CWarnerUssingDeviceNet() `CWarnerUssingDeviceNet ()`

Initializes a new instance of the `CWarnerUssingDeviceNet` class.

11.124.2.2 `~CWarnerUssingDeviceNet()` `virtual ~CWarnerUssingDeviceNet () [virtual]`

11.124.2.3 `!CWarnerUssingDeviceNet()` `!CWarnerUssingDeviceNet ()`


11.124.3 Property Documentation

11.124.3.1 **WarnerUssingFunction** `CWarnerUssingFunctionNet^ WarnerUssingFunction [get]`

11.125 CWarnerUssingFunctionNet Class Reference

`CWarnerUssingFunctionNet` is the class to control the Ussing device

Inheritance diagram for `CWarnerUssingFunctionNet`:

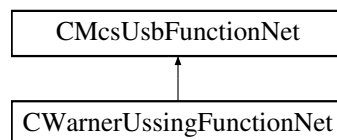

Public Member Functions

- `CWarnerUssingFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ pWarnerUssingFunctionPointerContainer)`
Initializes a new instance of the `CWarnerUssingFunctionNet` class.
- `CWarnerUssingFunctionNet (CMcsUsbNet^ mcsusb)`
- `virtual ~CWarnerUssingFunctionNet ()`
- `!CWarnerUssingFunctionNet ()`
- `int32_t GetChannelsCountOfChamber (int32_t ChamberId)`
gets number of channels in datastream from chamber amp with given index
- `int32_t GetNumberOfHardwareSlotsForChambers ()`
gets number of physical hardware slots for chambers amps
- `int32_t GetNumberOfAvailableChambers ()`
gets number of actually connected chamber amps
- `bool IsChamberAvailable (int32_t ChamberId)`
checks whether chamber amp is connected to slot
- `void SetPulse (int32_t ChamberId, UssingClampModeEnumNet StgMode, int32_t NumberOfRepetitions, array< int >^ Amplitudes, array< int >^ Durations, array< int >^ PulseMarker)`
defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes
- `void SetVoltageClampControllerParam_P (int32_t ChamberId, uint32_t P)`
sets P value of PID controller;
- `void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)`
sets I value of PID controller;

- void [SetVoltageClampControllerParam_D](#) (int32_t ChamberId, uint32_t D)
sets D value of PID controller;
- uint32_t [GetVoltageClampControllerParam_P](#) (int32_t ChamberId)
gets P value of PID controller;
- uint32_t [GetVoltageClampControllerParam_I](#) (int32_t ChamberId)
gets I value of PID controller;
- uint32_t [GetVoltageClampControllerParam_D](#) (int32_t ChamberId)
gets D value of PID controller;
- void [SetClampMode](#) (int32_t ChamberId, [UssingClampModeEnumNet](#) NewClampMode)
sets clamp mode (voltage, current or open clamp)
- [UssingClampModeEnumNet](#) [GetClampMode](#) (int32_t ChamberId)
gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)
- bool [IsInternalCalibrationFinished](#) (int32_t ChamberId)
when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available
- int32_t [GetU1Offset](#) (int32_t ChamberId)
- int32_t [GetU2Offset](#) (int32_t ChamberId)
- int32_t [GetU1Reference](#) (int32_t ChamberId)
- int32_t [GetU2Reference](#) (int32_t ChamberId)
- int32_t [GetDacZero](#) (int32_t ChamberId)
- void [SetHighCurrentMode](#) (int32_t ChamberId)
switch to high-current mode
- void [SetLowCurrentMode](#) (int32_t ChamberId)
switch to low-current mode
- bool [IsHighCurrentMode](#) (int32_t ChamberId)
- uint32_t [GetLowCurrentRange](#) (int32_t ChamberId)
query the range of the low current mode
- uint32_t [GetHighCurrentRange](#) (int32_t ChamberId)
query the range of the high current mode
- uint32_t [GetDacPampsPerDigitLowCurrentRange](#) (int32_t ChamberId)
get the resolution of the low current mode
- uint32_t [GetDacPampsPerDigitHighCurrentRange](#) (int32_t ChamberId)
get the resolution of the high current mode
- uint32_t [GetUnitsPerDigit](#) (int32_t ChamberId, int32_t ChannelId)
gets amps/volts per digit for specified chamber and channel
- int32_t [GetUnitExponent](#) (int32_t ChamberId, int32_t ChannelId)
gets the unit exponent for specified chamber and channel
- [UssingUnitEnumNet](#) [GetUnitName](#) (int32_t ChamberId, int32_t ChannelId)
gets the channel's unit name
- String ^ [GetUnitDescription](#) (int32_t ChamberId, int32_t ChannelId)
gets the description for the unit
- array< int > ^ [GetAvailableChambers](#) ()
returns array with (zero-based) ChamberIds of all available chambers
- int32_t [GetUptimeSeconds](#) (int32_t ChamberId)
gets the uptime in seconds
- void [SetIdleModeOffset](#) (int32_t ChamberId, [UssingClampModeEnumNet](#) ClampMode, int32_t NewIdleOffset)
sets the offset (voltage or current) that will be applied when clamping is DISABLED
- int32_t [GetIdleModeOffset](#) (int32_t ChamberId, [UssingClampModeEnumNet](#) ClampMode)
gets the offset (voltage or current) that will be applied when clamping is DISABLED
- void [SetEnablePulse](#) (int32_t ChamberId, [UssingClampModeEnumNet](#) ClampMode, bool Enable)
enable pulse of given chamber and mode (voltage/current clamp) of this chamber

- bool [IsPulseEnabled](#) (int32_t ChamberId, [UssingClampModeEnumNet](#) ClampMode)
returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED
- void [SetLiquidResistance](#) (int32_t ChamberId, int32_t NewLiquidResistance_Ohm)
sets the resistance of the liquid
- int32_t [GetLiquidResistance](#) (int32_t ChamberId)
gets the resistance of the liquid
- int32_t [GetComplianceVoltageThreshold](#) (int32_t ChamberId)
returns compliance voltage threshold in uV; when Uc is above, current source is overloaded
- bool [CompensateElectrodeOffset](#) (int32_t ChamberId)
blocking call to compensate electrode offset of one chamber; returns true when successful
- bool [WaitForChamber](#) (int32_t ChamberId)
blocking call that waits for chamber boot-up calibration to complete
- bool [WaitForAllChambers](#) ()
blocking call that waits for ALL chambers' boot-up calibration to complete

Additional Inherited Members

11.125.1 Detailed Description

[CWarnerUssingFunctionNet](#) is the class to control the Ussing device

11.125.2 Constructor & Destructor Documentation

11.125.2.1 CWarnerUssingFunctionNet() [1/2] [CWarnerUssingFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb,
[CMcsUsbFunctionPointerContainer](#)[^] pWarnerUssingFunctionPointerContainer)

Initializes a new instance of the [CWarnerUssingFunctionNet](#) class.

11.125.2.2 CWarnerUssingFunctionNet() [2/2] [CWarnerUssingFunctionNet](#) (
[CMcsUsbNet](#)[^] mcsusb)

11.125.2.3 ~CWarnerUssingFunctionNet() virtual [~CWarnerUssingFunctionNet](#) () [virtual]

11.125.2.4 "!CWarnerUssingFunctionNet() [!CWarnerUssingFunctionNet](#) ()

11.125.3 Member Function Documentation

11.125.3.1 CompensateElectrodeOffset() bool [CompensateElectrodeOffset](#) (
int32_t ChamberId)

blocking call to compensate electrode offset of one chamber; returns true when successful

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

true if compensation succeeded

11.125.3.2 GetAvailableChambers() `array<int> ^ GetAvailableChambers ()`

returns array with (zero-based) ChamberIds of all available chambers

11.125.3.3 GetChannelsCountOfChamber() `int32_t GetChannelsCountOfChamber (
int32_t ChamberId)`

gets number of channels in datastream from chamber amp with given index

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

return value of zero means that amp is not placed

11.125.3.4 GetClampMode() `UssingClampModeEnumNet GetClampMode (
int32_t ChamberId)`

gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

the current clamp mode

11.125.3.5 GetComplianceVoltageThreshold() `int32_t GetComplianceVoltageThreshold (`
`int32_t ChamberId)`

returns compliance voltage threshold in uV; when U_c is above, current source is overloaded

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

the compliance voltage threshold in uV

11.125.3.6 GetDacPampsPerDigitHighCurrentRange() `uint32_t GetDacPampsPerDigitHighCurrentRange`
`(`
`int32_t ChamberId)`

get the resolution of the high current mode

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

unit: pA/digit in high current mode

11.125.3.7 GetDacPampsPerDigitLowCurrentRange() `uint32_t GetDacPampsPerDigitLowCurrentRange (`
`int32_t ChamberId)`

get the resolution of the low current mode

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

pA/digit in low current mode

11.125.3.8 GetDacZero() `int32_t GetDacZero (`
`int32_t ChamberId)`

- diagnostic function only - ; gets real zero value of DAC in digits (0 -> neg. current; 32767 -> near zero; 65535 -> pos. current)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

the zero value of the DAC

11.125.3.9 GetHighCurrentRange() `uint32_t GetHighCurrentRange (`
`int32_t ChamberId)`

query the range of the high current mode

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

low current range in nA

11.125.3.10 GetIdleModeOffset() `int32_t GetIdleModeOffset (`
`int32_t ChamberId,`
`UssingClampModeEnumNet ClampMode)`

gets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation

Returns

unit: nA or uV

11.125.3.11 GetLiquidResistance() `int32_t GetLiquidResistance (`
`int32_t ChamberId)`

gets the resistance of the liquid

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

the liquid resistance in ohm

11.125.3.12 GetLowCurrentRange() `uint32_t GetLowCurrentRange (`
`int32_t ChamberId)`

query the range of the low current mode

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

low current range in nA

11.125.3.13 GetNumberOfAvailableChambers() `int32_t GetNumberOfAvailableChambers ()`

gets number of actually connected chamber amps

Returns

the number of actually connected chambers

11.125.3.14 GetNumberOfHardwareSlotsForChambers() `int32_t GetNumberOfHardwareSlotsForChambers ()`

gets number of physical hardware slots for chambers amps

Returns

the number of hardware chamber slots on the backplane

11.125.3.15 GetU1Offset() `int32_t GetU1Offset (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

U1 offset

11.125.3.16 GetU1Reference() `int32_t GetU1Reference (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

U1 reference

11.125.3.17 GetU2Offset() `int32_t GetU2Offset (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

U2 offset

11.125.3.18 GetU2Reference() `int32_t GetU2Reference (`
`int32_t ChamberId)`

- diagnostic function only -

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

U2 reference

11.125.3.19 GetUnitDescription() `String ^ GetUnitDescription (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets the description for the unit

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

the description of the unix

11.125.3.20 GetUnitExponent() `int32_t GetUnitExponent (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets the unit exponent for specified chamber and channel

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

example: return value -9 means that amps/volts per digit is in nano

11.125.3.21 GetUnitName() `UssingUnitEnumNet GetUnitName (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets the channel's unit name

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

the name of the unit

11.125.3.22 GetUnitsPerDigit() `uint32_t GetUnitsPerDigit (`
 `int32_t ChamberId,`
 `int32_t ChannelId)`

gets amps/volts per digit for specified chamber and channel

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ChannelId</i>	index of channel (zero-based)

Returns

amps/volts per digit

11.125.3.23 GetUptimeSeconds() `int32_t GetUptimeSeconds (`
 `int32_t ChamberId)`

gets the uptime in seconds

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

seconds since power-on

11.125.3.24 GetVoltageClampControllerParam_D() `uint32_t GetVoltageClampControllerParam_D (`
`int32_t ChamberId)`

gets D value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

the D value

11.125.3.25 GetVoltageClampControllerParam_I() `uint32_t GetVoltageClampControllerParam_I (`
`int32_t ChamberId)`

gets I value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

the I value

11.125.3.26 GetVoltageClampControllerParam_P() `uint32_t GetVoltageClampControllerParam_P (`
`int32_t ChamberId)`

gets P value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

Returns

the P value

11.125.3.27 IsChamberAvailable() `bool IsChamberAvailable (`
`int32_t ChamberId)`

checks whether chamber amp is connected to slot

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

true if the chamber is available

11.125.3.28 IsHighCurrentMode() `bool IsHighCurrentMode (`
`int32_t ChamberId)`

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

true if in high current mode

11.125.3.29 IsInternalCalibrationFinished() `bool IsInternalCalibrationFinished (`
`int32_t ChamberId)`

when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---

Returns

true if finished

11.125.3.30 IsPulseEnabled() `bool IsPulseEnabled (`
`int32_t ChamberId,`
`UssingClampModeEnumNet ClampMode)`

returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation

Returns

when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

11.125.3.31 SetClampMode() `void SetClampMode (`
 `int32_t ChamberId,`
 `UssingClampModeEnumNet NewClampMode)`

sets clamp mode (voltage, current or open clamp)

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>NewClampMode</i>	the clamp mode to use

11.125.3.32 SetEnablePulse() `void SetEnablePulse (`
 `int32_t ChamberId,`
 `UssingClampModeEnumNet ClampMode,`
 `bool Enable)`

enable pulse of given chamber and mode (voltage/current clamp) of this chamber

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation
<i>Enable</i>	when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

11.125.3.33 SetHighCurrentMode() `void SetHighCurrentMode (`
 `int32_t ChamberId)`

switch to high-current mode

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

11.125.3.34 SetIdleModeOffset() void SetIdleModeOffset (
 int32_t *ChamberId*,
 UssingClampModeEnumNet *ClampMode*,
 int32_t *NewIdleOffset*)

sets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>ClampMode</i>	voltage or current clamp stimulation
<i>NewIdleOffset</i>	unit: nA or uV

11.125.3.35 SetLiquidResistance() void SetLiquidResistance (
 int32_t *ChamberId*,
 int32_t *NewLiquidResistance_Ohm*)

sets the resistance of the liquid

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>NewLiquidResistance_Ohm</i>	the liquid resistance in ohm

11.125.3.36 SetLowCurrentMode() void SetLowCurrentMode (
 int32_t *ChamberId*)

switch to low-current mode

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
------------------------	---

11.125.3.37 SetPulse() void SetPulse (
 int32_t *ChamberId*,

```

UssingClampModeEnumNet StgMode,
int32_t NumberOfRepetitions,
array< int >^ Amplitudes,
array< int >^ Durations,
array< int >^ PulseMarker )

```

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based); send pattern to connected amp
<i>StgMode</i>	voltage or current clamp stimulation
<i>NumberOfRepetitions</i>	number of repetitions for pulse pattern (-1 for infinite; n means pattern is applied n+1 times)
<i>Amplitudes</i>	amplitude; unit in voltage clamp: uV; unit in current clamp: nA
<i>Durations</i>	duration in 100us; CAUTION: first element is applied only one; auto-loop back to second element after last one
<i>PulseMarker</i>	defines values on digital channel for each step (positive: digital channel "01", neg: "10", zero: "00")

11.125.3.38 SetVoltageClampControllerParam_D() void SetVoltageClampControllerParam_D (
int32_t ChamberId,
uint32_t D)

sets D value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>D</i>	useful range: 0..700

11.125.3.39 SetVoltageClampControllerParam_I() void SetVoltageClampControllerParam_I (
int32_t ChamberId,
uint32_t I)

sets I value of PID controller;

Parameters

<i>Chamber↔ Id</i>	index of hardware chamber slot (zero-based)
<i>I</i>	useful range: 80000..120000

11.125.3.40 SetVoltageClampControllerParam_P() void SetVoltageClampControllerParam_P (
 int32_t ChamberId,
 uint32_t P)

sets P value of PID controller;

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
<i>P</i>	useful value: 130000

11.125.3.41 WaitForAllChambers() bool WaitForAllChambers ()

blocking call that waits for ALL chambers' boot-up calibration to complete

Returns

returns false when at least one chamber's calibration fails (e.g. timeout...)

11.125.3.42 WaitForChamber() bool WaitForChamber (
 int32_t ChamberId)

blocking call that waits for chamber boot-up calibration to complete

Parameters

<i>ChamberId</i>	index of hardware chamber slot (zero-based)
------------------	---


Returns

returns false when calibration fails (e.g. timeout...)

11.126 CWarnerValveControllerDeviceNet Class Reference

[CWarnerValveControllerDeviceNet](#) is the class to access the Warner Valve Controller

Inheritance diagram for CWarnerValveControllerDeviceNet:

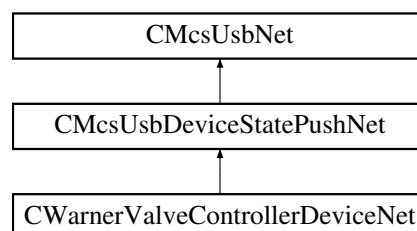

Public Member Functions

- delegate void [OnGetValveActive](#) (uint16_t valve, int valveActive)
- delegate void [OnGetValveManualState](#) (uint16_t valve, int32_t valveManualState)
- delegate void [OnGetValveManualGroup](#) (uint16_t valve, int32_t valveManualGroup)
- delegate void [OnGetValveMode](#) (uint16_t valve, [WvcValveModeEnumNet](#) ValveMode)
- delegate void [OnGetAnalogThresholdLow](#) (uint16_t valve, int32_t threshold)
- delegate void [OnGetAnalogThresholdHigh](#) (uint16_t valve, int32_t threshold)
- delegate void [OnGetDigitalPortDirection](#) (uint16_t port, [PortDirectionEnumNet](#) direction)
- delegate void [OnIsValveDigitalInInverted](#) (uint16_t valve, bool isInverted)
- delegate void [OnGetValveDigitalInPort](#) (uint16_t valve, uint32_t digitalInPort)
- delegate void [OnIsDigitalOutPortInverted](#) (uint16_t digitalOutPort, bool isInverted)
- delegate void [OnGetDigitalOutPortValve](#) (uint16_t digitalOutPort, uint32_t valve)
- delegate void [OnIsValveOpen](#) (uint16_t valve, bool valveOpen)
- delegate void [OnIsValveOpenInDigitalMode](#) (uint16_t valve, bool valveOpen)
- delegate void [OnIsValveOpenInAnalogMode](#) (uint16_t valve, bool valveOpen)
- delegate void [OnGetAnalogVoltage](#) (int32_t voltage)
- delegate void [OnTableEntryChanged](#) (uint16_t tableNumber)
- delegate void [OnGetTableNamebyIndex](#) (uint16_t tableNumber, String^ tableName)
- delegate void [OnGetActiveRunningTableNumber](#) (uint32_t tableNumber)
- delegate void [OnGetCurrentNumberOfValves](#) (int32_t numberOfValves)
- delegate void [OnGetValveBoardRevision](#) (uint32_t revision)
- delegate void [OnGetValveLedOn](#) (bool ledon)
- delegate void [OnGetDisplayMode](#) ([WvcDisplayModeEnumNet](#) DisplayMode)
- [CWarnerValveControllerDeviceNet](#) ()
Initializes a new instance of the [CWarnerValveControllerDeviceNet](#) class.
- virtual [~CWarnerValveControllerDeviceNet](#) ()
- [ICWarnerValveControllerDeviceNet](#) ()
- int [GetValveActive](#) (uint16_t valve)
Gets the valve active/inactive state
- void [SetValveActive](#) (uint16_t valve, int valveActive)
Sets the valve active/inactive state
- uint32_t [GetValvesActiveMap](#) ()
Gets the valves active/inactive states
- void [SetValvesActiveMap](#) (uint32_t valvesActive)
Sets the valve active/inactive state
- int32_t [GetValveManualState](#) (uint16_t valve)
Gets the valve manual on/off state
- void [SetValveManualState](#) (uint16_t valve, int32_t valveManualState)
Sets the valve manual on/off state
- uint32_t [GetValvesManualStateMap](#) ()
Gets the valves manual on/off states
- void [SetValvesManualStateMap](#) (uint32_t valveaManualState)
Sets the valve manual on/off state
- int32_t [GetValveManualGroup](#) (uint16_t valve)
Gets the valve manual group
- void [SetValveManualGroup](#) (uint16_t valve, int32_t valveManualGroup)
Sets the valve manual group
- [WvcValveModeEnumNet](#) [GetValveMode](#) (uint16_t valve)
Reads the valve mode
- void [SetValveMode](#) (uint16_t valve, [WvcValveModeEnumNet](#) ValveMode)
Writes the valve mode
- int32_t [GetAnalogThresholdLow](#) (uint16_t valve)

- Gets the lower threshold for the analog in port per valve*
- void [SetAnalogThresholdLow](#) (uint16_t valve, int32_t threshold)
- Sets the lower threshold for the analog in port per valve*
- int32_t [GetAnalogThresholdHigh](#) (uint16_t valve)
- Gets the upper threshold for the analog in port per valve*
- void [SetAnalogThresholdHigh](#) (uint16_t valve, int32_t threshold)
- Sets the upper threshold for the analog in port per valve*
- [PortDirectionEnumNet](#) [GetDigitalPortDirection](#) (uint16_t port)
- Gets the direction of a digital port*
- void [SetDigitalPortDirection](#) (uint16_t port, [PortDirectionEnumNet](#) direction)
- Sets the direction of a digital port*
- bool [IsValveDigitalInInverted](#) (uint16_t valve)
- Is digital in inverted*
- void [SetValveDigitalInInvert](#) (uint16_t valve, bool isInverted)
- Invert digital in*
- uint32_t [GetValveDigitalInPort](#) (uint16_t valve)
- Gets the number of the digital in port which is mapped to a valve*
- void [SetValveDigitalInPort](#) (uint16_t valve, uint32_t digitalInPort)
- Map a digital in port to a valve*
- bool [IsDigitalOutPortInverted](#) (uint16_t digitalOutPort)
- Gets the number of the valve which is mapped to a digital out port*
- void [SetDigitalOutPortInvert](#) (uint16_t digitalOutPort, bool isInverted)
- Map a valve to a digital out port*
- uint32_t [GetDigitalOutPortValve](#) (uint16_t digitalOutPort)
- Gets the number of the valve which is mapped to a digital out port*
- void [SetDigitalOutPortValve](#) (uint16_t digitalOutPort, uint32_t valve)
- Map a valve to a digital out port*
- void [SetDefault](#) ()
- Sets the settings of the valve controller to default*
- bool [IsValveOpen](#) (uint16_t valve)
- Is valve open*
- bool [IsValveOpenInDigitalMode](#) (uint16_t valve)
- True, if the valve would be open when the device is in digital mode*
- bool [IsValveOpenInAnalogMode](#) (uint16_t valve)
- True, if the valve would be open when the device is in analog mode*
- int32_t [GetAnalogVoltage](#) ()
- Reads the voltage on the analog in port*
- void [GetValveTableEntry](#) (uint16_t valve, uint16_t index, [System::Runtime::InteropServices::Out]uint32_t% duration, [System::Runtime::InteropServices::Out]bool% state)
- Read an entry from the valve protocol table*
- void [SetValveTableEntry](#) (uint16_t valve, uint16_t index, uint32_t duration, bool state)
- Write an entry to the valve protocol table*
- void [ClearValveTable](#) (uint16_t valve)
- Clear the valve protocol table*
- void [LoadValveTable](#) ()
- Load the current table from permanent memory*
- void [StoreValveTable](#) ()
- Store the current table in permanent memory*
- String ^ [GetTableNamebyIndex](#) (uint16_t tableNumber)
- Get the name of a protocol table*
- String ^ [GetTableName](#) ()

- Get the name of the current protocol table*
 - void [SetTableName](#) (String^ tableName)
- Set the name of the current protocol table*
 - uint32_t [GetActiveRunningTableNumber](#) ()
- Gets the number of the table that is active for running*
 - void [SetActiveRunningTableNumber](#) (uint32_t tableNumber)
- Sets the number of the table that is active for running*
 - uint32_t [GetCurrentEditTableNumber](#) ()
- Gets the number of the table that is current for editing*
 - void [SetCurrentEditTableNumber](#) (uint32_t tableNumber)
- Sets the number of the table that is current for editing*
 - void [ClearTableName](#) ()
- Clear the name of current protocol table*
 - void [SetTableStep](#) (uint16_t valve, int32_t steps)
- Skips the table protocol for a valve by steps*
 - void [SetTableStepAll](#) (int32_t steps)
- Skips the table protocol for all valves by steps*
 - int32_t [GetTotalNumberOfValves](#) ()
- Get the total number of valves in the system*
 - int32_t [GetTotalNumberOfDigitalPorts](#) ()
- Get the total number of digital ports in the system*
 - int32_t [GetTotalTableSize](#) ()
- Get the total table size in the system*
 - int32_t [GetTotalNumberOfTables](#) ()
- Get the total number of tables in the system*
 - int32_t [GetCurrentNumberOfValves](#) ()
- Get the current number of valves connected to the system*
 - uint32_t [GetValveBoardRevision](#) ()
- Gets the revision code of the valve board*
 - bool [GetValveLedOn](#) ()
- Gets the LED state of the valve board*
 - void [SetValveLedOn](#) (bool ledon)
- Gets the LED state of the valve board*
 - [WvcDisplayModeEnumNet](#) [GetDisplayMode](#) ()
- Reads the display mode*
 - void [SetDisplayMode](#) ([WvcDisplayModeEnumNet](#) DisplayMode, int32_t lockTimeMs)
- Writes the display mode*
 - String^ [GetValveBoardRevisionString](#) ()
- Gets the revision name of the valve board*
 - void [GetValveCurrent](#) ([System::Runtime::InteropServices::Out]int16_t% switch_current, [System::Runtime::InteropServices::Out]int16_t% hold_current)
- Gets the valve currents*
 - void [SetValveCurrent](#) (int16_t switch_current, int16_t hold_current)
- Sets the valve currents different from the default*

Events

- [OnGetValveActive](#)[^] [GetValveActiveEvent](#) [add, remove, raise]
Event fires when the valve state for the valve number has changed
- [OnGetValveManualState](#)[^] [GetValveManualStateEvent](#) [add, remove, raise]
Event fires when the manual valve state for the valve number has changed
- [OnGetValveManualGroup](#)[^] [GetValveManualGroupEvent](#) [add, remove, raise]
Event fires when the manual valve group for the valve number has changed
- [OnGetValveMode](#)[^] [GetValveModeEvent](#) [add, remove, raise]
Event fires when the valve mode for the valve number has changed
- [OnGetAnalogThresholdLow](#)[^] [GetAnalogThresholdLowEvent](#) [add, remove, raise]
Event fires when the threshold in mV for the valve number has changed
- [OnGetAnalogThresholdHigh](#)[^] [GetAnalogThresholdHighEvent](#) [add, remove, raise]
Event fires when the threshold in mV for the valve number has changed
- [OnGetDigitalPortDirection](#)[^] [GetDigitalPortDirectionEvent](#) [add, remove, raise]
Event fires when the direction for the port number has changed
- [OnIsValveDigitalInInverted](#)[^] [IsValveDigitalInInvertedEvent](#) [add, remove, raise]
Event fires when is inverted for the valve number has changed
- [OnGetValveDigitalInPort](#)[^] [GetValveDigitalInPortEvent](#) [add, remove, raise]
Event fires when the digital in port for the valve number has changed
- [OnIsDigitalOutPortInverted](#)[^] [IsDigitalOutPortInvertedEvent](#) [add, remove, raise]
Event fires when is inverted for the digital out port has changed
- [OnGetDigitalOutPortValve](#)[^] [GetDigitalOutPortValveEvent](#) [add, remove, raise]
Event fires when the valve number for the digital out port has changed
- [OnIsValveOpen](#)[^] [IsValveOpenEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnIsValveOpenInDigitalMode](#)[^] [IsValveOpenInDigitalModeEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnIsValveOpenInAnalogMode](#)[^] [IsValveOpenInAnalogModeEvent](#) [add, remove, raise]
Event fires when is open for the valve number has changed
- [OnGetAnalogVoltage](#)[^] [GetAnalogVoltageEvent](#) [add, remove, raise]
Event fires when the voltage in mV has changed
- [OnTableEntryChanged](#)[^] [TableEntryChangedEvent](#) [add, remove, raise]
Event fires when an entry of a table changed
- [OnGetTableNamebyIndex](#)[^] [GetTableNamebyIndexEvent](#) [add, remove, raise]
Event fires when the name of the table for the table number has changed
- [OnGetActiveRunningTableNumber](#)[^] [GetActiveRunningTableNumberEvent](#) [add, remove, raise]
Event fires when the table number has changed
- [OnGetCurrentNumberOfValves](#)[^] [GetCurrentNumberOfValvesEvent](#) [add, remove, raise]
Event fires when the number of valves has changed
- [OnGetValveBoardRevision](#)[^] [GetValveBoardRevisionEvent](#) [add, remove, raise]
Event fires when the revision code has changed
- [OnGetValveLedOn](#)[^] [GetValveLedOnEvent](#) [add, remove, raise]
Event fires when the LED state has changed
- [OnGetDisplayMode](#)[^] [GetDisplayModeEvent](#) [add, remove, raise]
Event fires when the display mode has changed

Additional Inherited Members

11.126.1 Detailed Description

[CWarnerValveControllerDeviceNet](#) is the class to access the Warner Valve Controller

11.126.2 Constructor & Destructor Documentation

11.126.2.1 CWarnerValveControllerDeviceNet() CWarnerValveControllerDeviceNet ()

Initializes a new instance of the CWarnerValveControllerDeviceNet class.

11.126.2.2 ~CWarnerValveControllerDeviceNet() virtual ~CWarnerValveControllerDeviceNet () [virtual]

11.126.2.3 !CWarnerValveControllerDeviceNet() !CWarnerValveControllerDeviceNet ()

11.126.3 Member Function Documentation

11.126.3.1 ClearTableName() void ClearTableName ()

Clear the name of current protocol table

11.126.3.2 ClearValveTable() void ClearValveTable (uint16_t valve)

Clear the valve protocol table

Parameters

<i>valve</i>	The valve number
--------------	------------------

11.126.3.3 GetActiveRunningTableNumber() uint32_t GetActiveRunningTableNumber ()

Gets the number of the table that is active for running

Returns

The table number

11.126.3.4 GetAnalogThresholdHigh() `int32_t GetAnalogThresholdHigh (uint16_t valve)`

Gets the upper threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The threshold in mV

11.126.3.5 GetAnalogThresholdLow() `int32_t GetAnalogThresholdLow (uint16_t valve)`

Gets the lower threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The threshold in mV

11.126.3.6 GetAnalogVoltage() `int32_t GetAnalogVoltage ()`

Reads the voltage on the analog in port

Returns

The voltage in mV

11.126.3.7 GetCurrentEditTableNumber() `uint32_t GetCurrentEditTableNumber ()`

Gets the number of the table that is current for editing

Returns

The table number

11.126.3.8 GetCurrentNumberOfValves() `int32_t GetCurrentNumberOfValves ()`

Get the current number of valves connected to the system

Returns

The number of valves

11.126.3.9 GetDigitalOutPortValve() `uint32_t GetDigitalOutPortValve (uint16_t digitalOutPort)`

Gets the number of the valve which is mapped to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
-----------------------	----------------------

Returns

The valve number

11.126.3.10 GetDigitalPortDirection() `PortDirectionEnumNet GetDigitalPortDirection (uint16_t port)`

Gets the direction of a digital port

Parameters

<i>port</i>	The port number
-------------	-----------------

Returns

the direction

11.126.3.11 GetDisplayMode() `WvcDisplayModeEnumNet GetDisplayMode ()`

Reads the display mode

Returns

The display mode

11.126.3.12 GetTableName() `String ^ GetTableName ()`

Get the name of the current protocol table

Returns

The name of the table

11.126.3.13 GetTableNamebyIndex() `String ^ GetTableNamebyIndex (
uint16_t tableNumber)`

Get the name of a protocol table

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

Returns

The name of the table

11.126.3.14 GetTotalNumberOfDigitalPorts() `int32_t GetTotalNumberOfDigitalPorts ()`

Get the total number of digital ports in the system

Returns

The number of digital ports

11.126.3.15 GetTotalNumberOfTables() `int32_t GetTotalNumberOfTables ()`

Get the total number of tables in the system

Returns

The number of tables

11.126.3.16 GetTotalNumberOfValves() `int32_t GetTotalNumberOfValves ()`

Get the total number of valves in the system

Returns

The number of valves

11.126.3.17 GetTotalTableSize() `int32_t GetTotalTableSize ()`

Get the total table size in the system

Returns

The table size

11.126.3.18 GetValveActive() `int GetValveActive (
uint16_t valve)`

Gets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve state

11.126.3.19 GetValveBoardRevision() `uint32_t GetValveBoardRevision ()`

Gets the revision code of the valve board

Returns

The revision code

11.126.3.20 GetValveBoardRevisionString() `String ^ GetValveBoardRevisionString ()`

Gets the revision name of the valve board

Returns

The revision name

11.126.3.21 GetValveCurrent() void GetValveCurrent (
 [System::Runtime::InteropServices::Out] int16_t% *switch_current*,
 [System::Runtime::InteropServices::Out] int16_t% *hold_current*)

Gets the valve currents

Parameters

<i>switch_current</i>	The switch current (in DAC units)
<i>hold_current</i>	The hold current (in DAC units)

11.126.3.22 GetValveDigitalInPort() uint32_t GetValveDigitalInPort (
 uint16_t *valve*)

Gets the number of the digital in port which is mapped to a valve

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The digital in port

11.126.3.23 GetValveLedOn() bool GetValveLedOn ()

Gets the LED state of the valve board

Returns

The LED state

11.126.3.24 GetValveManualGroup() int32_t GetValveManualGroup (
 uint16_t *valve*)

Gets the valve manual group

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The manual valve group

11.126.3.25 GetValveManualState() `int32_t GetValveManualState (`
`uint16_t valve)`

Gets the valve manual on/off state

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The manual valve state

11.126.3.26 GetValveMode() `WvcValveModeEnumNet GetValveMode (`
`uint16_t valve)`

Reads the valve mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

The valve mode

11.126.3.27 GetValvesActiveMap() `uint32_t GetValvesActiveMap ()`

Gets the valves active/inactive states

Returns

The valves states

11.126.3.28 GetValvesManualStateMap() `uint32_t GetValvesManualStateMap ()`

Gets the valves manual on/off states

Returns

The manual valves states

11.126.3.29 GetValveTableEntry() `void GetValveTableEntry (
 uint16_t valve,
 uint16_t index,
 [System::Runtime::InteropServices::Out] uint32_t% duration,
 [System::Runtime::InteropServices::Out] bool% state)`

Read an entry from the valve protocol table

Parameters

<i>valve</i>	The valve number
<i>index</i>	The index in the table
<i>duration</i>	the duration in ms
<i>state</i>	the state

11.126.3.30 IsDigitalOutPortInverted() `bool IsDigitalOutPortInverted (
 uint16_t digitalOutPort)`

Gets the number of the valve which is mapped to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
-----------------------	----------------------

Returns

is inverted

11.126.3.31 IsValveDigitalInInverted() `bool IsValveDigitalInInverted (
 uint16_t valve)`

Is digital in inverted

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is inverted

11.126.3.32 IsValveOpen() `bool IsValveOpen (`
 `uint16_t valve)`

Is valve open

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

11.126.3.33 IsValveOpenInAnalogMode() `bool IsValveOpenInAnalogMode (`
 `uint16_t valve)`

True, if the valve would be open when the device is in analog mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

11.126.3.34 IsValveOpenInDigitalMode() `bool IsValveOpenInDigitalMode (`
 `uint16_t valve)`

True, if the valve would be open when the device is in digital mode

Parameters

<i>valve</i>	The valve number
--------------	------------------

Returns

is open

11.126.3.35 LoadValveTable() `void LoadValveTable ()`

Load the current table from permanent memory

11.126.3.36 OnGetActiveRunningTableNumber() `delegate void OnGetActiveRunningTableNumber (uint32_t tableNumber)`**11.126.3.37 OnGetAnalogThresholdHigh()** `delegate void OnGetAnalogThresholdHigh (uint16_t valve, int32_t threshold)`**11.126.3.38 OnGetAnalogThresholdLow()** `delegate void OnGetAnalogThresholdLow (uint16_t valve, int32_t threshold)`**11.126.3.39 OnGetAnalogVoltage()** `delegate void OnGetAnalogVoltage (int32_t voltage)`**11.126.3.40 OnGetCurrentNumberOfValves()** `delegate void OnGetCurrentNumberOfValves (int32_t numberOfValves)`**11.126.3.41 OnGetDigitalOutPortValve()** `delegate void OnGetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)`**11.126.3.42 OnGetDigitalPortDirection()** `delegate void OnGetDigitalPortDirection (uint16_t port, PortDirectionEnumNet direction)`**11.126.3.43 OnGetDisplayMode()** `delegate void OnGetDisplayMode (WvcDisplayModeEnumNet DisplayMode)`

11.126.3.44 OnGetTableNamebyIndex() `delegate void OnGetTableNamebyIndex (`
 `uint16_t tableNumber,`
 `String^ tableName)`

11.126.3.45 OnGetValveActive() `delegate void OnGetValveActive (`
 `uint16_t valve,`
 `int valveActive)`

11.126.3.46 OnGetValveBoardRevision() `delegate void OnGetValveBoardRevision (`
 `uint32_t revision)`

11.126.3.47 OnGetValveDigitalInPort() `delegate void OnGetValveDigitalInPort (`
 `uint16_t valve,`
 `uint32_t digitalInPort)`

11.126.3.48 OnGetValveLedOn() `delegate void OnGetValveLedOn (`
 `bool ledon)`

11.126.3.49 OnGetValveManualGroup() `delegate void OnGetValveManualGroup (`
 `uint16_t valve,`
 `int32_t valveManualGroup)`

11.126.3.50 OnGetValveManualState() `delegate void OnGetValveManualState (`
 `uint16_t valve,`
 `int32_t valveManualState)`

11.126.3.51 OnGetValveMode() `delegate void OnGetValveMode (`
 `uint16_t valve,`
 `WvcValveModeEnumNet ValveMode)`

11.126.3.52 OnIsDigitalOutPortInverted() delegate void OnIsDigitalOutPortInverted (
 uint16_t digitalOutPort,
 bool isInverted)

11.126.3.53 OnIsValveDigitalInInverted() delegate void OnIsValveDigitalInInverted (
 uint16_t valve,
 bool isInverted)

11.126.3.54 OnIsValveOpen() delegate void OnIsValveOpen (
 uint16_t valve,
 bool valveOpen)

11.126.3.55 OnIsValveOpenInAnalogMode() delegate void OnIsValveOpenInAnalogMode (
 uint16_t valve,
 bool valveOpen)

11.126.3.56 OnIsValveOpenInDigitalMode() delegate void OnIsValveOpenInDigitalMode (
 uint16_t valve,
 bool valveOpen)

11.126.3.57 OnTableEntryChanged() delegate void OnTableEntryChanged (
 uint16_t tableNumber)

11.126.3.58 SetActiveRunningTableNumber() void SetActiveRunningTableNumber (
 uint32_t tableNumber)

Sets the number of the tanle that is active for running

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

11.126.3.59 SetAnalogThresholdHigh() void SetAnalogThresholdHigh (
 uint16_t valve,
 int32_t threshold)

Sets the upper threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
<i>threshold</i>	The threshold in mV

11.126.3.60 SetAnalogThresholdLow() `void SetAnalogThresholdLow (`
 `uint16_t valve,`
 `int32_t threshold)`

Sets the lower threshold for the analog in port per valve

Parameters

<i>valve</i>	The valve number
<i>threshold</i>	The threshold in mV

11.126.3.61 SetCurrentEditTableNumber() `void SetCurrentEditTableNumber (`
 `uint32_t tableNumber)`

Sets the number of the table that is current for editing

Parameters

<i>tableNumber</i>	The table number
--------------------	------------------

11.126.3.62 SetDefault() `void SetDefault ()`

Sets the settings of the valve controller to default

11.126.3.63 SetDigitalOutPortInvert() `void SetDigitalOutPortInvert (`
 `uint16_t digitalOutPort,`
 `bool isInverted)`

Map a valve to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
<i>isInverted</i>	True if digital out is to be inverted

11.126.3.64 SetDigitalOutPortValve() `void SetDigitalOutPortValve (`
 `uint16_t digitalOutPort,`
 `uint32_t valve)`

Map a valve to a digital out port

Parameters

<i>digitalOutPort</i>	The digital out port
<i>valve</i>	The valve number

11.126.3.65 SetDigitalPortDirection() `void SetDigitalPortDirection (`
 `uint16_t port,`
 `PortDirectionEnumNet direction)`

Sets the direction of a digital port

Parameters

<i>port</i>	The port number
<i>direction</i>	the direction

11.126.3.66 SetDisplayMode() `void SetDisplayMode (`
 `WvcDisplayModeEnumNet DisplayMode,`
 `int32_t lockTimeMs)`

Writes the display mode

Parameters

<i>DisplayMode</i>	The display mode
<i>lockTimeMs</i>	Locks the display for ms

11.126.3.67 SetTableName() `void SetTableName (`
 `String^ tableName)`

Set the name of the current protocol table

Parameters

<i>tableName</i>	The name of the table
------------------	-----------------------

11.126.3.68 SetTableStep() `void SetTableStep (`
 `uint16_t valve,`
 `int32_t steps)`

Skips the table protocol for a valve by steps

Parameters

<i>valve</i>	The valve number
<i>steps</i>	Number of steps

11.126.3.69 SetTableStepAll() `void SetTableStepAll (`
 `int32_t steps)`

Skips the table protocol for all valves by steps

Parameters

<i>steps</i>	Number of steps
--------------	-----------------

11.126.3.70 SetValveActive() `void SetValveActive (`
 `uint16_t valve,`
 `int valveActive)`

Sets the valve active/inactive state

Parameters

<i>valve</i>	The valve number
<i>valveActive</i>	The valve state

11.126.3.71 SetValveCurrent() `void SetValveCurrent (`
 `int16_t switch_current,`
 `int16_t hold_current)`

Sets the valve currents different from the default

Parameters

<i>switch_current</i>	The switch current (in DAC units); -1 sets the device default current
<i>hold_current</i>	The hold current (in DAC units); -1 sets the device default current

11.126.3.72 SetValveDigitalInInvert() `void SetValveDigitalInInvert (`
 `uint16_t valve,`
 `bool isInverted)`

Invert digital in

Parameters

<i>valve</i>	The valve number
<i>isInverted</i>	True if digital in is to be inverted

11.126.3.73 SetValveDigitalInPort() `void SetValveDigitalInPort (`
 `uint16_t valve,`
 `uint32_t digitalInPort)`

Map a digital in port to a valve

Parameters

<i>valve</i>	The valve number
<i>digitalInPort</i>	The digital in port

11.126.3.74 SetValveLedOn() `void SetValveLedOn (`
 `bool ledon)`

Gets the LED state of the valve board

Parameters

<i>ledon</i>	The LED state
--------------	---------------

11.126.3.75 SetValveManualGroup() `void SetValveManualGroup (`
 `uint16_t valve,`
 `int32_t valveManualGroup)`

Sets the valve manual group

Parameters

<i>valve</i>	The valve number
<i>valveManualGroup</i>	The manual valve group

11.126.3.76 SetValveManualState() `void SetValveManualState (`
 `uint16_t valve,`
 `int32_t valveManualState)`

Sets the valve manual on/off state

Parameters

<i>valve</i>	The valve number
<i>valveManualState</i>	The manual valve state

11.126.3.77 SetValveMode() `void SetValveMode (`
 `uint16_t valve,`
 `WvcValveModeEnumNet ValveMode)`

Writes the valve mode

Parameters

<i>valve</i>	The valve number
<i>ValveMode</i>	The valve mode

11.126.3.78 SetValvesActiveMap() `void SetValvesActiveMap (`
 `uint32_t valvesActive)`

Sets the valve active/inactive state

Parameters

<i>valvesActive</i>	The valves states
---------------------	-------------------

11.126.3.79 SetValvesManualStateMap() `void SetValvesManualStateMap (`
 `uint32_t valveaManualState)`

Sets the valve manual on/off state

Parameters

<i>valveaManualState</i>	The manual valves states
--------------------------	--------------------------

11.126.3.80 SetValveTableEntry() `void SetValveTableEntry (`
 `uint16_t valve,`
 `uint16_t index,`
 `uint32_t duration,`
 `bool state)`

Write an entry to the valve protocol table

Parameters

<i>valve</i>	The valve number
<i>index</i>	The index in the table
<i>duration</i>	the duration in ms
<i>state</i>	the state

11.126.3.81 StoreValveTable() `void StoreValveTable ()`

Store the current table in permanent memory

11.126.4 Event Documentation

11.126.4.1 GetActiveRunningTableNumberEvent `OnGetActiveRunningTableNumber^ GetActiveRunning↔`
`TableNumberEvent [add], [remove], [raise]`

Event fires when the table number has changed

11.126.4.2 GetAnalogThresholdHighEvent `OnGetAnalogThresholdHigh^ GetAnalogThresholdHighEvent`
`[add], [remove], [raise]`

Event fires when the threshold in mV for the valve number has changed

11.126.4.3 GetAnalogThresholdLowEvent `OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent`
`[add], [remove], [raise]`

Event fires when the threshold in mV for the valve number has changed

11.126.4.4 GetAnalogVoltageEvent `OnGetAnalogVoltage^ GetAnalogVoltageEvent [add], [remove], [raise]`

Event fires when the voltage in mV has changed

11.126.4.5 GetCurrentNumberOfValvesEvent `OnGetCurrentNumberOfValves^ GetCurrentNumberOfValvesEvent [add], [remove], [raise]`

Event fires when the number of valves has changed

11.126.4.6 GetDigitalOutPortValveEvent `OnGetDigitalOutPortValve^ GetDigitalOutPortValveEvent [add], [remove], [raise]`

Event fires when the valve number for the digital out port has changed

11.126.4.7 GetDigitalPortDirectionEvent `OnGetDigitalPortDirection^ GetDigitalPortDirectionEvent [add], [remove], [raise]`

Event fires when the direction for the port number has changed

11.126.4.8 GetDisplayModeEvent `OnGetDisplayMode^ GetDisplayModeEvent [add], [remove], [raise]`

Event fires when the display mode has changed

11.126.4.9 GetTableNamebyIndexEvent `OnGetTableNamebyIndex^ GetTableNamebyIndexEvent [add], [remove], [raise]`

Event fires when the name of the table for the table number has changed

11.126.4.10 GetValveActiveEvent `OnGetValveActive^ GetValveActiveEvent [add], [remove], [raise]`

Event fires when the valve state for the valve number has changed

11.126.4.11 GetValveBoardRevisionEvent [OnGetValveBoardRevision](#)[^] GetValveBoardRevisionEvent
[add], [remove], [raise]

Event fires when the revision code has changed

11.126.4.12 GetValveDigitalInPortEvent [OnGetValveDigitalInPort](#)[^] GetValveDigitalInPortEvent [add],
[remove], [raise]

Event fires when the digital in port for the valve number has changed

11.126.4.13 GetValveLedOnEvent [OnGetValveLedOn](#)[^] GetValveLedOnEvent [add], [remove], [raise]

Event fires when the LED state has changed

11.126.4.14 GetValveManualGroupEvent [OnGetValveManualGroup](#)[^] GetValveManualGroupEvent [add],
[remove], [raise]

Event fires when the manual valve group for the valve number has changed

11.126.4.15 GetValveManualStateEvent [OnGetValveManualState](#)[^] GetValveManualStateEvent [add],
[remove], [raise]

Event fires when the manual valve state for the valve number has changed

11.126.4.16 GetValveModeEvent [OnGetValveMode](#)[^] GetValveModeEvent [add], [remove], [raise]

Event fires when the valve mode for the valve number has changed

11.126.4.17 IsDigitalOutPortInvertedEvent [OnIsDigitalOutPortInverted](#)[^] IsDigitalOutPortInverted↔
Event [add], [remove], [raise]

Event fires when is inverted for the digital out port has changed

11.126.4.18 IsValveDigitalInInvertedEvent [OnIsValveDigitalInInverted](#)[^] IsValveDigitalInInverted↔
Event [add], [remove], [raise]

Event fires when is inverted for the valve number has changed

11.126.4.19 IsValveOpenEvent [OnIsValveOpen](#)[^] IsValveOpenEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.126.4.20 IsValveOpenInAnalogModeEvent [OnIsValveOpenInAnalogMode](#)[^] IsValveOpenInAnalogMode↔
Event [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.126.4.21 IsValveOpenInDigitalModeEvent [OnIsValveOpenInDigitalMode](#)[^] IsValveOpenInDigital↔
ModeEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed


11.126.4.22 TableEntryChangedEvent [OnTableEntryChanged](#)[^] TableEntryChangedEvent [add], [remove], [raise]

Event fires when an entry of a table changed

11.127 CWarnerValveControllerDeviceTesterFunctionNet Class Reference

[CWarnerValveControllerDeviceTesterFunctionNet](#) is the class to access the functions for the Warner Valve Controller Device Tester

Inheritance diagram for CWarnerValveControllerDeviceTesterFunctionNet:

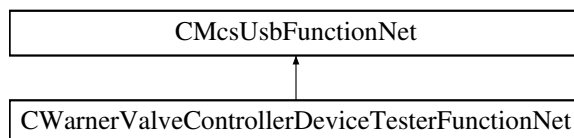

Public Member Functions

- [CWarnerValveControllerDeviceTesterFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] pWarnerValveControllerDeviceTesterFunctionPointerContainer)
Initializes a new instance of the [CWarnerValveControllerDeviceTesterFunctionNet](#) class.
- [CWarnerValveControllerDeviceTesterFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb)
- virtual [~CWarnerValveControllerDeviceTesterFunctionNet](#) ()
- [!CWarnerValveControllerDeviceTesterFunctionNet](#) ()
- void [SetADC](#) (uint32_t onoff)
Sets the ADC port of the tester
- uint32_t [GetSync](#) ()
Gets the output from the sync port
- void [SetTrigger](#) (uint32_t trigger)
Sets the input to the trigger port
- void [SetTriggerSyncDirection](#) (uint32_t direction)
Sets the direction of the trigger/sync test port
- uint32_t [GetIO](#) ()
Gets the output from the io ports
- void [SetIO](#) (uint32_t io)
Sets the input to the io ports
- void [SetIODirection](#) (int32_t direction)
Sets the direction of the IO test ports

Additional Inherited Members

11.127.1 Detailed Description

[CWarnerValveControllerDeviceTesterFunctionNet](#) is the class to access the functions for the Warner Valve Controller Device Tester

11.127.2 Constructor & Destructor Documentation

11.127.2.1 CWarnerValveControllerDeviceTesterFunctionNet() [1/2] [CWarnerValveControllerDeviceTesterFunctionNet](#)
 (
 [CMcsUsbNet](#)[^] mcsusb,
 [CMcsUsbFunctionPointerContainer](#)[^] pWarnerValveControllerDeviceTesterFunction↔
 PointerContainer)

Initializes a new instance of the [CWarnerValveControllerDeviceTesterFunctionNet](#) class.

11.127.2.2 CWarnerValveControllerDeviceTesterFunctionNet() [2/2] [CWarnerValveControllerDeviceTesterFunctionNet](#)
 (
 [CMcsUsbNet](#)[^] mcsusb)

11.127.2.3 `~CWarnerValveControllerDeviceTesterFunctionNet()` `virtual ~CWarnerValveControllerDeviceTesterFunctionNet()`
() [virtual]

11.127.2.4 `!CWarnerValveControllerDeviceTesterFunctionNet()` `!CWarnerValveControllerDeviceTesterFunctionNet()`
()

11.127.3 Member Function Documentation

11.127.3.1 `GetIO()` `uint32_t GetIO ()`

Gets the output from the io ports

Returns

The manual valves states

11.127.3.2 `GetSync()` `uint32_t GetSync ()`

Gets the output from the sync port

Returns

The sync state

11.127.3.3 `SetADC()` `void SetADC (`
`uint32_t onoff)`

Sets the ADC port of the tester

Parameters

<i>onoff</i>	The port state
--------------	----------------

11.127.3.4 `SetIO()` `void SetIO (`
`uint32_t io)`

Sets the input to the io ports

Parameters

<i>io</i>	The manual valves states
-----------	--------------------------

11.127.3.5 SetIODirection() `void SetIODirection (`
`int32_t direction)`

Sets the direction of the IO test ports

Parameters

<i>direction</i>	The 16bit direction map: 1=IN 0=OUT
------------------	-------------------------------------

11.127.3.6 SetTrigger() `void SetTrigger (`
`uint32_t trigger)`

Sets the input to the trigger port

Parameters

<i>trigger</i>	The trigger state
----------------	-------------------

11.127.3.7 SetTriggerSyncDirection() `void SetTriggerSyncDirection (`
`uint32_t direction)`


Sets the direction of the trigger/sync test port

Parameters

<i>direction</i>	The direction: 1=IN 0=OUT
------------------	---------------------------

11.128 CWClassicFunctionNet Class Reference

Inheritance diagram for CWClassicFunctionNet:

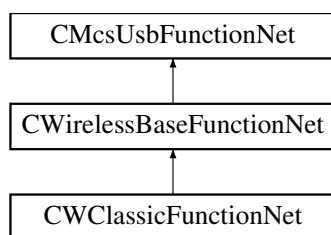

Public Member Functions

- [CWClassicFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb, [CMcsUsbFunctionPointerContainer](#)^ wClassicFuntion↔ PointerContainer)
- [CWClassicFunctionNet](#) ([CMcsUsbNet](#)^ mcsusb)
- [uint32_t ResetChannelmap](#) (unsigned int virtualDevice)
- [uint32_t SetChannelmap](#) (unsigned char position, unsigned char channel, unsigned int Device)
- void [SetHWSelectedChannels](#) (array< bool >^ channels, unsigned int Device)
- void [SetRFLostBehaviour](#) (uint8_t stoponfailure, unsigned int Device)
- void [SetHeadstageOnOff](#) (uint16_t onoff)
- USHORT [GetHeadstageOnOff](#) ()
- void [SetRFFrequencyHeadstage](#) (uint8_t receiver_nb, unsigned short frequency)
- unsigned short [GetRFFrequencyHeadstage](#) (uint8_t receiver_nb)
- void [SetRFFrequencyReceiver](#) (uint8_t receiver_nb, uint8_t configuration, unsigned short frequency)
- void [SetRFFrequencyReceiverEeprom](#) (uint8_t receiver_nb, uint8_t configuration, unsigned short frequency)
- unsigned short [GetRFFrequencyReceiver](#) (uint8_t receiver_nb, uint8_t configuration)
- void [SetSerialNumberHeadstage](#) (unsigned short number)
- unsigned short [GetSerialNumberHeadstage](#) ()
- void [SetSelectedHeadstage](#) (uint8_t number)
- [uint8_t GetSelectedHeadstage](#) ()
- void [ScanForHeadstages](#) ()
- [uint8_t GetScanHeadstagesResult](#) (int max_wait_for_ms)
- void [SetFilterParametersHeadstage](#) (unsigned short index, array< int >^ buffer)
- array< int > ^ [GetFilterParametersHeadstage](#) (unsigned short index)
- bool [GetHasRedLedHeadstage](#) ()
- void [SetHasChecksum](#) (unsigned int has, unsigned int Device)
- unsigned int [GetHasChecksum](#) (unsigned int Device)
- void [SetResetFilter](#) (unsigned int reset, unsigned int Device)
- unsigned int [GetResetFilter](#) (unsigned int Device)
- void [SetWPAType](#) (unsigned short type, unsigned int Device)
- unsigned short [GetWPAType](#) (unsigned int Device)
- void [SetWPADebugMode](#) (unsigned int mode, unsigned int Device)
- unsigned int [GetWPADebugMode](#) (unsigned int Device)
- void [SetRFPower](#) (unsigned short power)
- unsigned short [GetRFPower](#) ()
- unsigned int [GetRFConnectionStatus](#) ()

Additional Inherited Members

11.128.1 Constructor & Destructor Documentation

11.128.1.1 CWClassicFunctionNet() [1/2] `CWClassicFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ wClassicFuntionPointerContainer)`

11.128.1.2 CWClassicFunctionNet() [2/2] `CWClassicFunctionNet (CMcsUsbNet^ mcsusb)`

11.128.2 Member Function Documentation

11.128.2.1 GetFilterParametersHeadstage() `array<int> ^ GetFilterParametersHeadstage (unsigned short index)`

11.128.2.2 GetHasChecksum() `unsigned int GetHasChecksum (unsigned int Device)`

11.128.2.3 GetHasRedLedHeadstage() `bool GetHasRedLedHeadstage ()`

11.128.2.4 GetHeadstageOnOff() `USHORT GetHeadstageOnOff ()`

11.128.2.5 GetResetFilter() `unsigned int GetResetFilter (unsigned int Device)`

11.128.2.6 GetRFConnectionStatus() `unsigned int GetRFConnectionStatus ()`

11.128.2.7 GetRFFrequencyHeadstage() `unsigned short GetRFFrequencyHeadstage (uint8_t receiver_nb)`

11.128.2.8 GetRFFrequencyReceiver() unsigned short GetRFFrequencyReceiver (
 uint8_t *receiver_nb*,
 uint8_t *configuration*)

11.128.2.9 GetRFPower() unsigned short GetRFPower ()

11.128.2.10 GetScanHeadstagesResult() uint8_t GetScanHeadstagesResult (
 int *max_wait_for_ms*)

11.128.2.11 GetSelectedHeadstage() uint8_t GetSelectedHeadstage ()

11.128.2.12 GetSerialNumberHeadstage() unsigned short GetSerialNumberHeadstage ()

11.128.2.13 GetWPADebugMode() unsigned int GetWPADebugMode (
 unsigned int *Device*)

11.128.2.14 GetWPAType() unsigned short GetWPAType (
 unsigned int *Device*)

11.128.2.15 ResetChannelmap() uint32_t ResetChannelmap (
 unsigned int *virtualDevice*)

11.128.2.16 ScanForHeadstages() void ScanForHeadstages ()

11.128.2.17 SetChannelmap() uint32_t SetChannelmap (
 unsigned char *position*,
 unsigned char *channel*,
 unsigned int *Device*)

11.128.2.18 SetFilterParametersHeadstage() void SetFilterParametersHeadstage (
 unsigned short *index*,
 array< int >^ *buffer*)

11.128.2.19 SetHasChecksum() void SetHasChecksum (
 unsigned int *has*,
 unsigned int *Device*)

11.128.2.20 SetHeadstageOnOff() void SetHeadstageOnOff (
 uint16_t *onoff*)

11.128.2.21 SetHWSelectedChannels() void SetHWSelectedChannels (
 array< bool >^ *channels*,
 unsigned int *Device*)

11.128.2.22 SetResetFilter() void SetResetFilter (
 unsigned int *reset*,
 unsigned int *Device*)

11.128.2.23 SetRFFrequencyHeadstage() void SetRFFrequencyHeadstage (
 uint8_t *receiver_nb*,
 unsigned short *frequency*)

11.128.2.24 SetRFFrequencyReceiver() void SetRFFrequencyReceiver (
 uint8_t *receiver_nb*,
 uint8_t *configuration*,
 unsigned short *frequency*)

11.128.2.25 SetRFFrequencyReceiverEeprom() void SetRFFrequencyReceiverEeprom (
 uint8_t *receiver_nb*,
 uint8_t *configuration*,
 unsigned short *frequency*)

11.128.2.26 SetRFLostBehaviour() void SetRFLostBehaviour (
 uint8_t *stoponfailure*,
 unsigned int *Device*)

11.128.2.27 SetRFPower() void SetRFPower (
 unsigned short *power*)

11.128.2.28 SetSelectedHeadstage() void SetSelectedHeadstage (
 uint8_t *number*)


11.128.2.29 SetSerialNumberHeadstage() void SetSerialNumberHeadstage (
 unsigned short *number*)

11.128.2.30 SetWPADebugMode() void SetWPADebugMode (
 unsigned int *mode*,
 unsigned int *Device*)

11.128.2.31 SetWPAType() void SetWPAType (
 unsigned short *type*,
 unsigned int *Device*)

11.129 CWirelessBaseFunctionNet Class Reference

Inheritance diagram for CWirelessBaseFunctionNet:

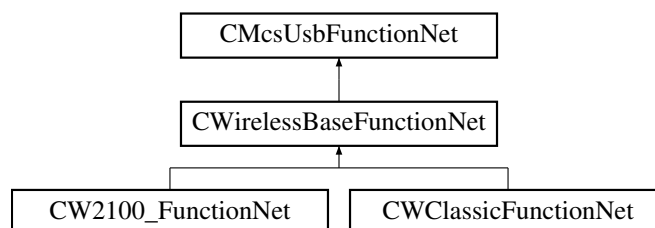

Public Member Functions

- [CWirelessBaseFunctionNet](#) ([CMcsUsbNet](#)[^] mcsusb, [CMcsUsbFunctionPointerContainer](#)[^] mcsusbfunction)

Static Public Member Functions

- static String ^ [CreateWirelessHeadstageSerialNumberString](#) (unsigned short ID)

Additional Inherited Members

11.129.1 Constructor & Destructor Documentation

11.129.1.1 CWirelessBaseFunctionNet() [CWirelessBaseFunctionNet](#) (
[CMcsUsbNet](#)^ *mcsusb*,
[CMcsUsbFunctionPointerContainer](#)^ *mcsusbfunction*)

11.129.2 Member Function Documentation

11.129.2.1 CreateWirelessHeadstageSerialNumberString() static String ^ [CreateWirelessHeadstageSerialNumberString](#) (
 unsigned short *ID*) [static]

11.130 DeviceldNet Struct Reference

Device Id.

Public Member Functions

- [DeviceldNet](#) ()
- [DeviceldNet](#) ([VendorIdEnumNet](#) vendor, [ProductIdEnumNet](#) product, int bcd, [McsBusTypeEnumNet](#) bustype)
- [DeviceldNet](#) ([DeviceldNet](#)% deviceId)
- [DeviceldNet](#) operator= ([DeviceldNet](#)% deviceId)

Public Attributes

- [VendorIdEnumNet](#) IdVendor
- [ProductIdEnumNet](#) IdProduct
- int [BcdDevice](#)
- [McsBusTypeEnumNet](#) BusType

11.130.1 Detailed Description

Device Id.

11.130.2 Constructor & Destructor Documentation

11.130.2.1 DeviceIdNet() [1/3] `DeviceIdNet ()`

11.130.2.2 DeviceIdNet() [2/3] `DeviceIdNet (`
 `VendorIdEnumNet vendor,`
 `ProductIdEnumNet product,`
 `int bcd,`
 `McsBusTypeEnumNet bustype)`

11.130.2.3 DeviceIdNet() [3/3] `DeviceIdNet (`
 `DeviceIdNet% deviceId)`

11.130.3 Member Function Documentation

11.130.3.1 operator=() `DeviceIdNet operator= (`
 `DeviceIdNet% deviceId)`

11.130.4 Member Data Documentation

11.130.4.1 BcdDevice `int BcdDevice`

11.130.4.2 BusType `McsBusTypeEnumNet BusType`

11.130.4.3 IdProduct `ProductIdEnumNet IdProduct`

11.130.4.4 IdVendor `VendorIdEnumNet IdVendor`

11.131 DigitalSource< digitalsourceenum > Class Template Reference

Public Member Functions

- [DigitalSource](#) ()
- [DigitalSource](#) (digitalsourceenum source)
- int [MaxBitNumber](#) ()
- int [MaxBitNumber](#) (digitalsourceenum [Source](#))

Static Public Member Functions

- static int [MaxBitNumberStatic](#) (digitalsourceenum [Source](#))
- static int [size](#) ()

Properties

- digitalsourceenum [Source](#) [get, set]

11.131.1 Constructor & Destructor Documentation

11.131.1.1 DigitalSource() [1/2] [DigitalSource](#) ()

11.131.1.2 DigitalSource() [2/2] [DigitalSource](#) (
digitalsourceenum *source*)

11.131.2 Member Function Documentation

11.131.2.1 MaxBitNumber() [1/2] int [MaxBitNumber](#) ()

11.131.2.2 MaxBitNumber() [2/2] int [MaxBitNumber](#) (
digitalsourceenum *Source*)

11.131.2.3 MaxBitNumberStatic() static int [MaxBitNumberStatic](#) (
digitalsourceenum *Source*) [static]

11.131.2.4 size() `static int size () [static]`

11.131.3 Property Documentation

11.131.3.1 Source `digitalsourceenum Source [get], [set]`

11.132 DigitalSourceGeneral Class Reference

Public Member Functions

- [DigitalSourceGeneral](#) (Type[^] type)
- [DigitalSourceGeneral](#) (Type[^] type, int [Source](#))
- int [MaxBitNumber](#) ()
- int [MaxBitNumber](#) (int [Source](#))

Static Public Member Functions

- static int [MaxBitNumber](#) (Type[^] type, int [Source](#))
- static int [size](#) (Type[^] type)

Properties

- int [Source](#) [get, set]

11.132.1 Constructor & Destructor Documentation

11.132.1.1 DigitalSourceGeneral() [1/2] [DigitalSourceGeneral](#) (
Type[^] type)

11.132.1.2 DigitalSourceGeneral() [2/2] [DigitalSourceGeneral](#) (
Type[^] type,
int *Source*)

11.132.2 Member Function Documentation

11.132.2.1 MaxBitNumber() [1/3] `int MaxBitNumber ()`

11.132.2.2 MaxBitNumber() [2/3] `int MaxBitNumber (`
`int Source)`

11.132.2.3 MaxBitNumber() [3/3] `static int MaxBitNumber (`
`Type^ type,`
`int Source) [static]`

11.132.2.4 size() `static int size (`
`Type^ type) [static]`

11.132.3 Property Documentation

11.132.3.1 Source `int Source [get], [set]`

11.133 DriverVersionNet Class Reference

Class gives firmware versions of the device's firmware destinations.

Public Member Functions

- [DriverVersionNet](#) ()
Constructor.
- [~DriverVersionNet](#) ()
Destructor.
- unsigned int [GetStatus](#) ([CFirmwareDestinationNet](#) dest)
Get status of firmware destination.
- unsigned int [GetStatus](#) (unsigned int index)
Get status of firmware destination.
- unsigned int [GetVersionInt](#) ([CFirmwareDestinationNet](#) dest)
Get the version number of firmware destination (major in high word, minor in low word)
- unsigned int [GetVersionInt](#) (unsigned int index)
Get the version number of firmware destination (major in high word, minor in low word)
- unsigned int [GetMajor](#) ([CFirmwareDestinationNet](#) dest)
Get the major version number of firmware destination.
- unsigned int [GetMajor](#) (unsigned int index)
Get the major version number of firmware destination.

- unsigned int [GetMinor](#) ([CFirmwareDestinationNet](#) dest)
Get the minor version number of firmware destination.
- unsigned int [GetMinor](#) (unsigned int index)
Get the minor version number of firmware destination.
- unsigned int [GetNumEntries](#) ()
Get the number of available firmware destinations.
- String ^ [GetVersionString](#) ([CFirmwareDestinationNet](#) dest)
Get the version as a string in the format Major.Minor.
- String ^ [GetVersionString](#) (unsigned int index)
Get the version as a string in the format Major.Minor.
- [CFirmwareDestinationNet](#) [GetDestinationCode](#) (unsigned int index)
Get CFirmwareDestinationNet.
- String ^ [GetDestinationName](#) ([CFirmwareDestinationNet](#) dest)
Get firmware destination name.
- String ^ [GetDestinationName](#) (unsigned int index)
Get firmware destination name.
- String ^ [GetSerialNumber](#) ([CFirmwareDestinationNet](#) dest)
Get the serial number of the destination, when no serial number if found, return an empty string.
- String ^ [GetSerialNumber](#) (unsigned int index)
Get the serial number of the destination, when no serial number if found, return an empty string.

Static Public Member Functions

- static String ^ [DriverVersionNet::FormatVersion](#) (unsigned int v)

11.133.1 Detailed Description

Class gives firmware versions of the device's firmware destinations.

11.133.2 Constructor & Destructor Documentation

11.133.2.1 [DriverVersionNet\(\)](#) [DriverVersionNet](#) ()

Constructor.

11.133.2.2 [~DriverVersionNet\(\)](#) [~DriverVersionNet](#) ()

Destructor.

11.133.3 Member Function Documentation

11.133.3.1 DriverVersionNet::FormatVersion() `static String ^ DriverVersionNet::FormatVersion (unsigned int v) [static]`

11.133.3.2 GetDestinationCode() `CFirmwareDestinationNet GetDestinationCode (unsigned int index)`

Get CFirmwareDestinationNet.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.3 GetDestinationName() [1/2] `String ^ GetDestinationName (CFirmwareDestinationNet dest)`

Get firmware destination name.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.4 GetDestinationName() [2/2] `String ^ GetDestinationName (unsigned int index)`

Get firmware destination name.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.5 GetMajor() [1/2] `unsigned int GetMajor (CFirmwareDestinationNet dest)`

Get the major version number of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.6 GetMajor() [2/2] `unsigned int GetMajor (unsigned int index)`

Get the major version number of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.7 GetMinor() [1/2] unsigned int GetMinor (
 CFirmwareDestinationNet dest)

Get the minor version number of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.8 GetMinor() [2/2] unsigned int GetMinor (
 unsigned int index)

Get the minor version number of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.9 GetNumEntries() unsigned int GetNumEntries ()

Get the number of available firmware destinations.

11.133.3.10 GetSerialNumber() [1/2] String ^ GetSerialNumber (
 CFirmwareDestinationNet dest)

Get the serial number of the destination, when no serial number if found, return an empty string.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.11 GetSerialNumber() [2/2] String ^ GetSerialNumber (
 unsigned int index)

Get the serial number of the destination, when no serial number if found, return an empty string.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.12 GetStatus() [1/2] unsigned int GetStatus (
CFirmwareDestinationNet *dest*)

Get status of firmware destination.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.13 GetStatus() [2/2] unsigned int GetStatus (
unsigned int *index*)

Get status of firmware destination.

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.14 GetVersionInt() [1/2] unsigned int GetVersionInt (
CFirmwareDestinationNet *dest*)

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.15 GetVersionInt() [2/2] unsigned int GetVersionInt (
unsigned int *index*)

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

<i>index</i>	by index of firmware destination
--------------	----------------------------------

11.133.3.16 GetVersionString() [1/2] `String ^ GetVersionString (CFirmwareDestinationNet dest)`

Get the version as a string in the format Major.Minor.

Parameters

<i>dest</i>	by CFirmwareDestinationNet
-------------	----------------------------

11.133.3.17 GetVersionString() [2/2] `String ^ GetVersionString (unsigned int index)`

Get the version as a string in the format Major.Minor.

Parameters

<i>index</i>	by index of firmware
--------------	----------------------

11.134 FirmwareDestinationNames Class Reference

Static Public Attributes

- static String ^ **DSP** = gcnew String("DSP")
- static String ^ **USB** = gcnew String("USB")
- static String ^ **MCU1** = gcnew String("MCU1")
- static String ^ **Bootstrap** = gcnew String("Bootstrap")
- static String ^ **MCSBUS1** = gcnew String("McsBus1")
- static String ^ **MCSBUS2** = gcnew String("McsBus2")
- static String ^ **MCSBUS3** = gcnew String("McsBus3")
- static String ^ **MCSBUS4** = gcnew String("McsBus4")
- static String ^ **MCSBUS5** = gcnew String("McsBus5")
- static String ^ **MCSBUS6** = gcnew String("McsBus6")
- static String ^ **MCSBUS7** = gcnew String("McsBus7")
- static String ^ **MCSBUS8** = gcnew String("McsBus8")
- static String ^ **MCSBUS9** = gcnew String("McsBus9")
- static String ^ **MCSBUS10** = gcnew String("McsBus10")
- static String ^ **MCSBUS11** = gcnew String("McsBus11")
- static String ^ **MCSBUS12** = gcnew String("McsBus12")
- static String ^ **MCSBUS13** = gcnew String("McsBus13")
- static String ^ **BUS1_MCSBUS1** = gcnew String("Bus1McsBus1")
- static String ^ **BUS1_MCSBUS2** = gcnew String("Bus1McsBus2")
- static String ^ **PIC** = gcnew String("PIC")
- static String ^ **PIC2** = gcnew String("PIC2")
- static String ^ **PIC3** = gcnew String("PIC3")
- static String ^ **PIC4** = gcnew String("PIC4")

- static String ^ Altera = gcnnew String("Altera")
- static String ^ FPGA2 = gcnnew String("FPGA2")
- static String ^ FPGA3 = gcnnew String("FPGA3")
- static String ^ FPGA4 = gcnnew String("FPGA4")
- static String ^ FPGA5 = gcnnew String("FPGA5")
- static String ^ FPGA6 = gcnnew String("FPGA6")

11.134.1 Member Data Documentation

11.134.1.1 Altera String ^ Altera = gcnnew String("Altera") [static]

11.134.1.2 Bootstrap String ^ Bootstrap = gcnnew String("Bootstrap") [static]

11.134.1.3 BUS1_MCSBUS1 String ^ BUS1_MCSBUS1 = gcnnew String("Bus1McsBus1") [static]

11.134.1.4 BUS1_MCSBUS2 String ^ BUS1_MCSBUS2 = gcnnew String("Bus1McsBus2") [static]

11.134.1.5 DSP String ^ DSP = gcnnew String("DSP") [static]

11.134.1.6 FPGA2 String ^ FPGA2 = gcnnew String("FPGA2") [static]

11.134.1.7 FPGA3 String ^ FPGA3 = gcnnew String("FPGA3") [static]

11.134.1.8 FPGA4 String ^ FPGA4 = gcnnew String("FPGA4") [static]

11.134.1.9 FPGA5 String ^ FPGA5 = gcnnew String("FPGA5") [static]

11.134.1.10 FPGA6 `String ^ FPGA6 = gcnew String("FPGA6") [static]`

11.134.1.11 MCSBUS1 `String ^ MCSBUS1 = gcnew String("McsBus1") [static]`

11.134.1.12 MCSBUS10 `String ^ MCSBUS10 = gcnew String("McsBus10") [static]`

11.134.1.13 MCSBUS11 `String ^ MCSBUS11 = gcnew String("McsBus11") [static]`

11.134.1.14 MCSBUS12 `String ^ MCSBUS12 = gcnew String("McsBus12") [static]`

11.134.1.15 MCSBUS13 `String ^ MCSBUS13 = gcnew String("McsBus13") [static]`

11.134.1.16 MCSBUS2 `String ^ MCSBUS2 = gcnew String("McsBus2") [static]`

11.134.1.17 MCSBUS3 `String ^ MCSBUS3 = gcnew String("McsBus3") [static]`

11.134.1.18 MCSBUS4 `String ^ MCSBUS4 = gcnew String("McsBus4") [static]`

11.134.1.19 MCSBUS5 `String ^ MCSBUS5 = gcnew String("McsBus5") [static]`

11.134.1.20 MCSBUS6 `String ^ MCSBUS6 = gcnew String("McsBus6") [static]`

11.134.1.21 MCSBUS7 `String ^ MCSBUS7 = gcnew String("McsBus7") [static]`

11.134.1.22 MCSBUS8 `String ^ MCSBUS8 = gcnew String("McsBus8") [static]`

11.134.1.23 MCSBUS9 `String ^ MCSBUS9 = gcnew String("McsBus9") [static]`

11.134.1.24 MCU1 `String ^ MCU1 = gcnew String("MCU1") [static]`

11.134.1.25 PIC `String ^ PIC = gcnew String("PIC") [static]`

11.134.1.26 PIC2 `String ^ PIC2 = gcnew String("PIC2") [static]`


11.134.1.27 PIC3 `String ^ PIC3 = gcnew String("PIC3") [static]`

11.134.1.28 PIC4 `String ^ PIC4 = gcnew String("PIC4") [static]`

11.134.1.29 USB `String ^ USB = gcnew String("USB") [static]`

11.135 HeadStageIDType Class Reference

Inheritance diagram for HeadStageIDType:

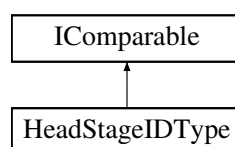

Public Types

- enum class [HeadstageTypeEnum](#) {
[Unknown](#) ,
[MeasuringOnly](#) ,
[OpticalStimulation](#) ,
[ElectricalStimulation](#) }

Public Member Functions

- [HeadStageIDType](#) (unsigned int entry, [CW2100_FunctionNet](#)^ device)
- virtual System::String ^ [ToString](#) () override
- virtual bool [Equals](#) (Object^ obj) override
- virtual Int32 [CompareTo](#) (Object^ obj)

Properties

- bool [Valid](#) [get]
- unsigned int [Entry](#) [get]
- unsigned short [ID](#) [get]
- System::String ^ [SN](#) [get]
- unsigned int [TypeValue](#) [get]
- System::String ^ [Type](#) [get]
- [HeadstageTypeEnum](#) [HeadstageType](#) [get]
- System::String ^ [UserDefinedName](#) [get]
- int [NumberOfAnalogChannels](#) [get]
- int [NumberOfStimulationChannels](#) [get]
- [W2100_StimulusParametersNet](#)^ [StimulusParameters](#) [get]
- bool [HasIMU](#) [get]
- bool [W16IsW14](#) [get]
- bool [HasOptoCurrentMeasurement](#) [get]

11.135.1 Member Enumeration Documentation

11.135.1.1 HeadstageTypeEnum enum [HeadstageTypeEnum](#) [strong]

Enumerator

Unknown	
MeasuringOnly	
OpticalStimulation	
ElectricalStimulation	

11.135.2 Constructor & Destructor Documentation

11.135.2.1 HeadStageIDType() `HeadStageIDType (`
 `unsigned int entry,`
 `CW2100_FunctionNet^ device)`

11.135.3 Member Function Documentation

11.135.3.1 CompareTo() `virtual Int32 CompareTo (`
 `Object^ obj) [virtual]`

11.135.3.2 Equals() `virtual bool Equals (`
 `Object^ obj) [override], [virtual]`

11.135.3.3 ToString() `virtual System::String ^ ToString () [override], [virtual]`

11.135.4 Property Documentation

11.135.4.1 Entry `unsigned int Entry [get]`

11.135.4.2 HasIMU `bool HasIMU [get]`

11.135.4.3 HasOptoCurrentMeasurement `bool HasOptoCurrentMeasurement [get]`

11.135.4.4 HeadstageType `HeadstageTypeEnum HeadstageType [get]`

11.135.4.5 ID `unsigned short ID [get]`

11.135.4.6 NumberOfAnalogChannels int NumberOfAnalogChannels [get]

11.135.4.7 NumberOfStimulationChannels int NumberOfStimulationChannels [get]

11.135.4.8 SN System:: String^ SN [get]

11.135.4.9 StimulusParameters W2100_StimulusParametersNet^ StimulusParameters [get]

11.135.4.10 Type System:: String^ Type [get]

11.135.4.11 TypeValue unsigned int TypeValue [get]

11.135.4.12 UserDefinedName System:: String^ UserDefinedName [get]

11.135.4.13 Valid bool Valid [get]

11.135.4.14 W16IsW14 bool W16IsW14 [get]

11.136 HeadstageIDTypeObject Class Reference

Public Member Functions

- [HeadstageIDTypeObject](#) ([HeadStageIDType](#)^ idType)
- virtual String ^ [ToString](#) () override
- virtual bool [Equals](#) (Object^ obj) override
- virtual int [GetHashCode](#) () override

Public Attributes

- [HeadStageIDType](#) ^ [_IdType](#)
- String ^ [_AdditionalText](#)

Properties

- [HeadStageIDType](#) ^ [IdType](#) [get]
- String ^ [AdditionalText](#) [get, set]

11.136.1 Constructor & Destructor Documentation

11.136.1.1 HeadstageIDTypeObject() [HeadstageIDTypeObject](#) (
 [HeadStageIDType](#) ^ *idType*)

11.136.2 Member Function Documentation

11.136.2.1 Equals() virtual bool Equals (
 Object ^ *obj*) [override], [virtual]

11.136.2.2 GetHashCode() virtual int GetHashCode () [override], [virtual]

11.136.2.3 ToString() virtual String ^ ToString () [override], [virtual]

11.136.3 Member Data Documentation

11.136.3.1 _AdditionalText String ^ [_AdditionalText](#)

11.136.3.2 _IdType [HeadStageIDType](#) ^ [_IdType](#)

11.136.4 Property Documentation

11.136.4.1 AdditionalText `String^ AdditionalText [get], [set]`

11.136.4.2 IdType `HeadStageIDType^ IdType [get]`

11.137 HeadStageIDTypeState Class Reference

Properties

- unsigned int [State](#) [get]
- [HeadStageIDType^](#) [IdType](#) [get]
- bool [ControlState](#) [get]
- bool [DataState](#) [get]

11.137.1 Property Documentation

11.137.1.1 ControlState `bool ControlState [get]`

11.137.1.2 DataState `bool DataState [get]`

11.137.1.3 IdType `HeadStageIDType^ IdType [get]`

11.137.1.4 State `unsigned int State [get]`

11.138 mkfilterNet Class Reference

Static Public Member Functions

- static int [mkfilter](#) (String[^] *filtertype*, double *value*, String[^] *passtype*, int *order*, double *alpha1*, double *alpha2*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Amplification*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*)
- static int [mkfilter_MCS_k](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Amplification*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *coeffs*)
- static int [mkfilter_MCS_k](#) (int *SamplesPerSecond*, double *R1*, double *R2*, double *C*, double *Correction*, [System::Runtime::InteropServices::Out] array< double >[^]% *coeffs*)
- static void [mkfilter_coef_in_one_set](#) (int *n*, [System::Runtime::InteropServices::In] array< double >[^] *xcoeffs*, [System::Runtime::InteropServices::In] array< double >[^] *ycoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *out_coeffs*)
- static void [mkfilter_scale_coef_in_one_set](#) (int *n*, double *scale*, [System::Runtime::InteropServices::In] array< double >[^] *xcoeffs*, [System::Runtime::InteropServices::In] array< double >[^] *ycoeffs*, [System::Runtime::InteropServices::Out] array< double >[^]% *out_coeffs*)
- static void [mkfilter_normalize_coeffs_short](#) (short *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< short >[^]% *out_coeffs*)
- static void [mkfilter_normalize_coeffs_int](#) (int *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< int >[^]% *out_coeffs*)
- static void [mkfilter_normalize_scale_coeffs_int](#) (int *maxvalue*, [System::Runtime::InteropServices::In] array< double >[^] *coeffs*, [System::Runtime::InteropServices::Out] array< int >[^]% *out_coeffs*)
- static double [mkfilter_highpass_coeff](#) (int *SamplesPerSecond*, double *Frequency*)
- static double [mkfilter_highpass_k](#) (int *SamplesPerSecond*, double *Frequency*)
- static double [mkfilter_highpass_frequency_from_coeff](#) (int *SamplesPerSecond*, double *coeff*)
- static double [mkfilter_highpass_frequency_from_k](#) (int *SamplesPerSecond*, double *k*)

11.138.1 Member Function Documentation

11.138.1.1 [mkfilter\(\)](#) static int [mkfilter](#) (
String[^] *filtertype*,
double *value*,
String[^] *passtype*,
int *order*,
double *alpha1*,
double *alpha2*,
[System::Runtime::InteropServices::Out] array< double >[^]% *xcoeffs*,
[System::Runtime::InteropServices::Out] array< double >[^]% *ycoeffs*) [static]

11.138.1.2 mkfilter_coef_in_one_set() static void mkfilter_coef_in_one_set (
 int *n*,
 [System::Runtime::InteropServices::In] array< double >^ *xcoeffs*,
 [System::Runtime::InteropServices::In] array< double >^ *ycoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *out_coeffs*) [static]

11.138.1.3 mkfilter_highpass_coeff() static double mkfilter_highpass_coeff (
 int *SamplesPerSecond*,
 double *Frequency*) [static]

11.138.1.4 mkfilter_highpass_frequency_from_coeff() static double mkfilter_highpass_frequency_from↵
 from_coeff (
 int *SamplesPerSecond*,
 double *coeff*) [static]

11.138.1.5 mkfilter_highpass_frequency_from_k() static double mkfilter_highpass_frequency_from↵
 _k (
 int *SamplesPerSecond*,
 double *k*) [static]

11.138.1.6 mkfilter_highpass_k() static double mkfilter_highpass_k (
 int *SamplesPerSecond*,
 double *Frequency*) [static]

11.138.1.7 mkfilter_MCS() [1/2] static int mkfilter_MCS (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Amplification*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *xcoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *ycoeffs*) [static]

11.138.1.8 mkfilter_MCS() [2/2] static int mkfilter_MCS (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *xcoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *ycoeffs*) [static]

- 11.138.1.9 mkfilter_MCS_k()** [1/2] static int mkfilter_MCS_k (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Amplification*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *coeffs*) [static]
- 11.138.1.10 mkfilter_MCS_k()** [2/2] static int mkfilter_MCS_k (
 int *SamplesPerSecond*,
 double *R1*,
 double *R2*,
 double *C*,
 double *Correction*,
 [System::Runtime::InteropServices::Out] array< double >^% *coeffs*) [static]
- 11.138.1.11 mkfilter_normalize_coeffs_int()** static void mkfilter_normalize_coeffs_int (
 int *maxvalue*,
 [System::Runtime::InteropServices::In] array< double >^ *coeffs*,
 [System::Runtime::InteropServices::Out] array< int >^% *out_coeffs*) [static]
- 11.138.1.12 mkfilter_normalize_coeffs_short()** static void mkfilter_normalize_coeffs_short (
 short *maxvalue*,
 [System::Runtime::InteropServices::In] array< double >^ *coeffs*,
 [System::Runtime::InteropServices::Out] array< short >^% *out_coeffs*) [static]
- 11.138.1.13 mkfilter_normalize_scale_coeffs_int()** static void mkfilter_normalize_scale_coeffs_int (
 int *maxvalue*,
 [System::Runtime::InteropServices::In] array< double >^ *coeffs*,
 [System::Runtime::InteropServices::Out] array< int >^% *out_coeffs*) [static]
- 11.138.1.14 mkfilter_scale_coef_in_one_set()** static void mkfilter_scale_coef_in_one_set (
 int *n*,
 double *scale*,
 [System::Runtime::InteropServices::In] array< double >^ *xcoeffs*,
 [System::Runtime::InteropServices::In] array< double >^ *ycoeffs*,
 [System::Runtime::InteropServices::Out] array< double >^% *out_coeffs*) [static]

11.139 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference

Public Member Functions

- void [SetParameter](#) (unsigned short command, unsigned short index, unsigned int value)
- void [SetParameter](#) (unsigned short command, unsigned short index, unsigned int value1, unsigned int value2)
- void [SetUserParameter](#) (unsigned short index, unsigned int value)
Stores persistently 32 bit integer values on RoboMain
- void [SetUserParameter](#) (unsigned short index, int value)
Stores persistently 32 bit integer values on RoboMain
- void [GetParameter](#) (unsigned short command, unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value)
- void [GetParameter](#) (unsigned short command, unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value1, [System::Runtime::InteropServices::Out]unsigned int% value2)
- void [GetUserParameter](#) (unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value)
Reads 32 bit integer values stored persistently on RoboMain
- void [GetUserParameter](#) (unsigned short index, [System::Runtime::InteropServices::Out]int% value)
Reads 32 bit integer values stored persistently on RoboMain
- void [FindReferencePhase0](#) (unsigned char busaddress, char axes)
- void [FindReferencePhase0](#) (unsigned char busaddress, char axes, int timeout)
- unsigned char [HasRef](#) (unsigned char busaddress, char axes)
- void [SetHWRevision](#) (unsigned int revision)
- unsigned int [GetHWRevision](#) ()
- void [SetHWConfig](#) (unsigned int config)
- unsigned int [GetHWConfig](#) ()
- void [SetMinPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMinPressureWaitTime](#) ()
- void [SetMinPressure](#) (unsigned int pressure)
- unsigned int [GetMinPressure](#) ()
- void [SetMaxPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMaxPressureWaitTime](#) ()
- void [SetMinNoPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMinNoPressureWaitTime](#) ()
- void [SetMaxNoPressure](#) (unsigned int pressure)
- unsigned int [GetMaxNoPressure](#) ()
- void [SetMaxNoPressureWaitTime](#) (unsigned int t)
- unsigned int [GetMaxNoPressureWaitTime](#) ()
- void [SetSearchReferenceMethod](#) (unsigned char busaddress, char axes, unsigned int method)
- unsigned int [GetSearchReferenceMethod](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceOffsetPos](#) (unsigned char busaddress, char axes, int offsetpos)
- int [GetSearchReferenceOffsetPos](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFastSpeed](#) (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short [GetSearchReferenceFastSpeed](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFastAccel](#) (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short [GetSearchReferenceFastAccel](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFineSpeed](#) (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short [GetSearchReferenceFineSpeed](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceFineAccel](#) (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short [GetSearchReferenceFineAccel](#) (unsigned char busaddress, char axes)
- void [SetSearchReferenceMoveOut](#) (unsigned char busaddress, char axes, int move)
- int [GetSearchReferenceMoveOut](#) (unsigned char busaddress, char axes)
- void [SetAxisConfig](#) (unsigned char busaddress, char axes, unsigned int config)
- unsigned int [GetAxisConfig](#) (unsigned char busaddress, char axes)
- void [GetPhases](#) (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out] unsigned short% phase0, [System::Runtime::InteropServices::Out] unsigned short% lastphase)

11.139.1 Member Function Documentation

11.139.1.1 FindReferencePhase0() [1/2] `void FindReferencePhase0 (`
 `unsigned char busaddress,`
 `char axes)`

11.139.1.2 FindReferencePhase0() [2/2] `void FindReferencePhase0 (`
 `unsigned char busaddress,`
 `char axes,`
 `int timeout)`

11.139.1.3 GetAxisConfig() `unsigned int GetAxisConfig (`
 `unsigned char busaddress,`
 `char axes)`

11.139.1.4 GetHWConfig() `unsigned int GetHWConfig ()`

11.139.1.5 GetHWRevision() `unsigned int GetHWRevision ()`

11.139.1.6 GetMaxNoPressure() `unsigned int GetMaxNoPressure ()`

11.139.1.7 GetMaxNoPressureWaitTime() `unsigned int GetMaxNoPressureWaitTime ()`

11.139.1.8 GetMaxPressureWaitTime() `unsigned int GetMaxPressureWaitTime ()`

11.139.1.9 GetMinNoPressureWaitTime() `unsigned int GetMinNoPressureWaitTime ()`

11.139.1.10 GetMinPressure() unsigned int GetMinPressure ()

11.139.1.11 GetMinPressureWaitTime() unsigned int GetMinPressureWaitTime ()

11.139.1.12 GetParameter() [1/2] void GetParameter (
 unsigned short *command*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned int% *value*)

11.139.1.13 GetParameter() [2/2] void GetParameter (
 unsigned short *command*,
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned int% *value1*,
 [System::Runtime::InteropServices::Out] unsigned int% *value2*)

11.139.1.14 GetPhases() void GetPhases (
 unsigned char *busaddress*,
 char *axes*,
 [System::Runtime::InteropServices::Out] unsigned short% *phase0*,
 [System::Runtime::InteropServices::Out] unsigned short% *lastphase*)

11.139.1.15 GetSearchReferenceFastAccel() unsigned short GetSearchReferenceFastAccel (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.16 GetSearchReferenceFastSpeed() unsigned short GetSearchReferenceFastSpeed (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.17 GetSearchReferenceFineAccel() unsigned short GetSearchReferenceFineAccel (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.18 GetSearchReferenceFineSpeed() unsigned short GetSearchReferenceFineSpeed (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.19 GetSearchReferenceMethod() unsigned int GetSearchReferenceMethod (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.20 GetSearchReferenceMoveOut() int GetSearchReferenceMoveOut (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.21 GetSearchReferenceOffsetPos() int GetSearchReferenceOffsetPos (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.22 GetUserParameter() [1/2] void GetUserParameter (
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] int% *value*)

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data buffer

11.139.1.23 GetUserParameter() [2/2] void GetUserParameter (
 unsigned short *index*,
 [System::Runtime::InteropServices::Out] unsigned int% *value*)

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data buffer

11.139.1.24 HasRef() unsigned char HasRef (
 unsigned char *busaddress*,
 char *axes*)

11.139.1.25 SetAxisConfig() void SetAxisConfig (
 unsigned char *busaddress*,
 char *axes*,
 unsigned int *config*)

11.139.1.26 SetHWConfig() void SetHWConfig (
 unsigned int *config*)

11.139.1.27 SetHWRevision() void SetHWRevision (
 unsigned int *revision*)

11.139.1.28 SetMaxNoPressure() void SetMaxNoPressure (
 unsigned int *pressure*)

11.139.1.29 SetMaxNoPressureWaitTime() void SetMaxNoPressureWaitTime (
 unsigned int *t*)

11.139.1.30 SetMaxPressureWaitTime() void SetMaxPressureWaitTime (
 unsigned int *t*)

11.139.1.31 SetMinNoPressureWaitTime() void SetMinNoPressureWaitTime (
 unsigned int *t*)

11.139.1.32 SetMinPressure() void SetMinPressure (
 unsigned int *pressure*)

11.139.1.33 SetMinPressureWaitTime() void SetMinPressureWaitTime (
 unsigned int *t*)

11.139.1.34 SetParameter() [1/2] void SetParameter (
 unsigned short *command*,
 unsigned short *index*,
 unsigned int *value*)

11.139.1.35 SetParameter() [2/2] void SetParameter (
 unsigned short *command*,
 unsigned short *index*,
 unsigned int *value1*,
 unsigned int *value2*)

11.139.1.36 SetSearchReferenceFastAccel() void SetSearchReferenceFastAccel (
 unsigned char *busaddress*,
 char *axes*,
 unsigned short *accel*)

11.139.1.37 SetSearchReferenceFastSpeed() void SetSearchReferenceFastSpeed (
 unsigned char *busaddress*,
 char *axes*,
 unsigned short *speed*)

11.139.1.38 SetSearchReferenceFineAccel() void SetSearchReferenceFineAccel (
 unsigned char *busaddress*,
 char *axes*,
 unsigned short *accel*)

11.139.1.39 SetSearchReferenceFineSpeed() void SetSearchReferenceFineSpeed (
 unsigned char *busaddress*,
 char *axes*,
 unsigned short *speed*)

11.139.1.40 SetSearchReferenceMethod() `void SetSearchReferenceMethod (`
 `unsigned char busaddress,`
 `char axes,`
 `unsigned int method)`

11.139.1.41 SetSearchReferenceMoveOut() `void SetSearchReferenceMoveOut (`
 `unsigned char busaddress,`
 `char axes,`
 `int move)`

11.139.1.42 SetSearchReferenceOffsetPos() `void SetSearchReferenceOffsetPos (`
 `unsigned char busaddress,`
 `char axes,`
 `int offsetpos)`

11.139.1.43 SetUserParameter() [1/2] `void SetUserParameter (`
 `unsigned short index,`
 `int value)`

Stores *persistently* 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data to be stored

11.139.1.44 SetUserParameter() [2/2] `void SetUserParameter (`
 `unsigned short index,`
 `unsigned int value)`

Stores *persistently* 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

<i>index</i>	address offset of parameter; range: 0..15
<i>value</i>	data to be stored

11.140 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference

Public Member Functions

- void [FindReferencePhase0XY](#) ()
- void [FindReferencePhase0XY](#) (int timeout)

11.140.1 Member Function Documentation

11.140.1.1 [FindReferencePhase0XY\(\)](#) [1/2] `void FindReferencePhase0XY ()`

11.140.1.2 [FindReferencePhase0XY\(\)](#) [2/2] `void FindReferencePhase0XY (
int timeout)`

11.141 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference

Public Member Functions

- [s_FilterAttributesNet](#) (s_FilterAttributes attrib)
- s_FilterAttributes [ToCpp](#) ()

Public Attributes

- uint32_t [PreCommaB](#)
- uint32_t [PostCommaB](#)
- uint32_t [CommaPositionB](#)
- uint32_t [PreCommaA](#)
- uint32_t [PostCommaA](#)
- uint32_t [CommaPositionA](#)

11.141.1 Constructor & Destructor Documentation

11.141.1.1 [s_FilterAttributesNet\(\)](#) `s_FilterAttributesNet (
s_FilterAttributes attrib)`

11.141.2 Member Function Documentation

11.141.2.1 ToCpp() `s_FilterAttributes ToCpp ()`

11.141.3 Member Data Documentation

11.141.3.1 CommaPositionA `uint32_t CommaPositionA`

11.141.3.2 CommaPositionB `uint32_t CommaPositionB`

11.141.3.3 PostCommaA `uint32_t PostCommaA`

11.141.3.4 PostCommaB `uint32_t PostCommaB`

11.141.3.5 PreCommaA `uint32_t PreCommaA`

11.141.3.6 PreCommaB `uint32_t PreCommaB`

11.142 CMeaAudioFunctionNet::s_setaudionet Struct Reference

Public Attributes

- int [channel](#)
- int [amplification](#)

11.142.1 Member Data Documentation

11.142.1.1 amplification `int amplification`

11.142.1.2 channel `int channel`

11.143 CStimulusFunctionNet::SidebandData Class Reference

Public Member Functions

- [SidebandData](#) ()
- [~SidebandData](#) ()
Destructor: called by Dispose()
- [!SidebandData](#) ()
Finalizer: called by GC before collecting

Properties

- `array< int32_t >^` [Sideband](#) [get]
- `array< uint64_t >^` [Duration](#) [get]

11.143.1 Constructor & Destructor Documentation

11.143.1.1 SidebandData() [SidebandData](#) ()

11.143.1.2 ~SidebandData() [~SidebandData](#) ()

Destructor: called by Dispose()

11.143.1.3 "!SidebandData() [!SidebandData](#) ()

Finalizer: called by GC before collecting

11.143.2 Property Documentation

11.143.2.1 Duration `array< uint64_t >^` [Duration](#) [get]

11.143.2.2 Sideband `array< int32_t >^` [Sideband](#) [get]

11.144 StgStatusNet Class Reference

Static Public Member Functions

- static [StgStatusNet](#) ^ [FromIntPtr](#) (IntPtr stgstatus)
- static [StgStatusNet](#) ^ [FromPtr](#) (stgstatus_t *stgstatus)

Public Attributes

- array< [Stg200xTriggerStatusEnumNet](#) > ^ [TiggerStatus](#)
- array< uint32_t > ^ [ListOfChangedTriggers](#)

11.144.1 Member Function Documentation

11.144.1.1 FromIntPtr() static [StgStatusNet](#) ^ FromIntPtr (IntPtr stgstatus) [static]

11.144.1.2 FromPtr() static [StgStatusNet](#) ^ FromPtr (stgstatus_t * stgstatus) [static]

11.144.2 Member Data Documentation

11.144.2.1 ListOfChangedTriggers array<uint32_t> ^ ListOfChangedTriggers

11.144.2.2 TiggerStatus array<[Stg200xTriggerStatusEnumNet](#)> ^ TiggerStatus

11.145 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference

Public Member Functions

- [StimulusDeviceDataAndUnrolledData](#) ()
- [~StimulusDeviceDataAndUnrolledData](#) ()
Destructor: called by Dispose()
- [!StimulusDeviceDataAndUnrolledData](#) ()
Finalizer: called by GC before collecting

Properties

- array< uint8_t >^ [DeviceData](#) [get]
- int [DeviceDataLength](#) [get]
- array< int32_t >^ [UnrolledAmplitude](#) [get]
- array< uint32_t >^ [UnrolledSync](#) [get]
- array< uint64_t >^ [UnrolledDuration](#) [get]

11.145.1 Constructor & Destructor Documentation

11.145.1.1 [StimulusDeviceDataAndUnrolledData\(\)](#) [StimulusDeviceDataAndUnrolledData](#) ()

11.145.1.2 [~StimulusDeviceDataAndUnrolledData\(\)](#) [~StimulusDeviceDataAndUnrolledData](#) ()

Destructor: called by Dispose()

11.145.1.3 ["!StimulusDeviceDataAndUnrolledData\(\)](#) [!StimulusDeviceDataAndUnrolledData](#) ()

Finalizer: called by GC before collecting

11.145.2 Property Documentation

11.145.2.1 [DeviceData](#) array< uint8_t >^ [DeviceData](#) [get]

11.145.2.2 [DeviceDataLength](#) int [DeviceDataLength](#) [get]

11.145.2.3 [UnrolledAmplitude](#) array< int32_t >^ [UnrolledAmplitude](#) [get]

11.145.2.4 [UnrolledDuration](#) array< uint64_t >^ [UnrolledDuration](#) [get]

11.145.2.5 UnrolledSync `array< uint32_t>^ UnrolledSync [get]`

11.146 usbSetupPacket_t Class Reference

Public Attributes

- `uint8_t` [bmRequestType](#)
- `uint8_t` [bRequest](#)
- `uint16_t` [wValue](#)
- `uint16_t` [wIndex](#)
- `uint16_t` [wLength](#)

11.146.1 Member Data Documentation

11.146.1.1 bmRequestType `uint8_t bmRequestType`

11.146.1.2 bRequest `uint8_t bRequest`

11.146.1.3 wIndex `uint16_t wIndex`

11.146.1.4 wLength `uint16_t wLength`

11.146.1.5 wValue `uint16_t wValue`

11.147 W2100_StimulusParametersNet Struct Reference

Public Attributes

- `int` [DACResolution](#)
- `int` [TimeResolutionInNanoSeconds](#)
- `int` [VoltageRangeInMicroVolt](#)
- `int` [VoltageResolutionInMicroVolt](#)
- `int` [CurrentRangeInNanoAmp](#)
- `int` [CurrentResolutionInNanoAmp](#)

11.147.1 Member Data Documentation

11.147.1.1 CurrentRangeInNanoAmp `int CurrentRangeInNanoAmp`

11.147.1.2 CurrentResolutionInNanoAmp `int CurrentResolutionInNanoAmp`

11.147.1.3 DACResolution `int DACResolution`

11.147.1.4 TimeResolutionInNanoSeconds `int TimeResolutionInNanoSeconds`

11.147.1.5 VoltageRangeInMicroVolt `int VoltageRangeInMicroVolt`

11.147.1.6 VoltageResolutionInMicroVolt `int VoltageResolutionInMicroVolt`

Index

- !CDacCalibrationFunctionNet
 - CDacCalibrationFunctionNet, [113](#)
- !CDigOutStimulatorFunctionNet
 - CDigOutStimulatorFunctionNet, [123](#)
- !CEXternDTesterDeviceNet
 - CEXternDTesterDeviceNet, [128](#)
- !CGrapheneFunctionNet
 - CGrapheneFunctionNet, [168](#)
- !CInterfaceboard2FunctionNet
 - CInterfaceboard2FunctionNet, [185](#)
- !CInterfaceboardFunctionNet
 - CInterfaceboardFunctionNet, [187](#)
- !CLIH3DeviceNet
 - CLIH3DeviceNet, [190](#)
- !CMEA2100x256FunctionNet
 - CMEA2100x256FunctionNet, [330](#)
- !CMcsUsbFunctionNet
 - CMcsUsbFunctionNet, [299](#)
- !CMcsUsbListNet
 - CMcsUsbListNet, [306](#)
- !CMcsUsbNet
 - CMcsUsbNet, [312](#)
- !CMeFunctionNet
 - CMeFunctionNet, [363](#)
- !CMeaCleanDeviceNet
 - CMeaCleanDeviceNet, [335](#)
- !CMeaCoatDeviceNet
 - CMeaCoatDeviceNet, [339](#)
- !CMultiBatteryChargerDeviceNet
 - CMultiBatteryChargerDeviceNet, [365](#)
- !CMultiwellCallbackFunctionNet
 - CMultiwellCallbackFunctionNet, [373](#)
- !CMultiwellDeviceNet
 - CMultiwellDeviceNet, [376](#)
- !CMultiwellOptoStimFunctionNet
 - CMultiwellOptoStimFunctionNet, [383](#)
- !CPPCFunctionNet
 - CPPCFunctionNet, [416](#)
- !CPedoterDeviceNet
 - CPedoterDeviceNet, [397](#)
- !CPositionIIDeviceNet
 - CPositionIIDeviceNet, [403](#)
- !CPositionImpDeviceNet
 - CPositionImpDeviceNet, [411](#)
- !CProgramPressureCurveNet
 - CProgramPressureCurveNet, [430](#)
- !CPulseGeneratorFunctionNet
 - CPulseGeneratorFunctionNet, [431](#)
- !CRFFFunctionNet
 - CRFFFunctionNet, [439](#)
- !CSCUFunctionNet
 - CSCUFunctionNet, [495](#)
- !CTEERFunctionNet
 - CTEERFunctionNet, [585](#)
- !CUsbDeviceConfigurationFunctionNet, [595](#)
- !CWarnerUssingDeviceNet
 - CWarnerUssingDeviceNet, [613](#)
- !CWarnerUssingFunctionNet
 - CWarnerUssingFunctionNet, [615](#)
- !CWarnerValveControllerDeviceNet
 - CWarnerValveControllerDeviceNet, [634](#)
- !CWarnerValveControllerDeviceTesterFunctionNet
 - CWarnerValveControllerDeviceTesterFunctionNet, [656](#)
- !SidebandData
 - CStimulusFunctionNet::SidebandData, [694](#)
- !StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- _AdditionalText
 - HeadstageIDTypeObject, [680](#)
- _IdType
 - HeadstageIDTypeObject, [680](#)
- ~CCMOSMeaDeviceNet
 - CCMOSMeaDeviceNet, [108](#)
- ~CChannelTestDeviceNet
 - CChannelTestDeviceNet, [96](#)
- ~CCreateFilterNet
 - CCreateFilterNet, [111](#)
- ~CDacCalibrationFunctionNet
 - CDacCalibrationFunctionNet, [113](#)
- ~CDigOutStimulatorFunctionNet
 - CDigOutStimulatorFunctionNet, [123](#)
- ~CEXternDTesterDeviceNet
 - CEXternDTesterDeviceNet, [128](#)
- ~CFilterCoefficientsNet
 - CFilterCoefficientsNet, [130](#)
- ~CFilterPropertyNet
 - CFilterPropertyNet, [135](#)
- ~CFluidControlDeviceNet
 - CFluidControlDeviceNet, [138](#)
- ~CGenericDevelopDeviceNet
 - CGenericDevelopDeviceNet, [152](#)
- ~CGilsonDeviceNet
 - CGilsonDeviceNet, [164](#)
- ~CGrapheneASICDeviceNet
 - CGrapheneASICDeviceNet, [165](#)
- ~CGrapheneFunctionNet
 - CGrapheneFunctionNet, [168](#)
- ~CInterfaceboard2FunctionNet
 - CInterfaceboard2FunctionNet, [185](#)
- ~CInterfaceboardFunctionNet
 - CInterfaceboardFunctionNet, [187](#)
- ~CLIH3DeviceNet
 - CLIH3DeviceNet, [190](#)
- ~CMEA2100x256FunctionNet
 - CMEA2100x256FunctionNet, [330](#)
- ~CMcsBusNet
 - CMcsBusNet, [236](#)

- ~CMcsBus_AxisParametersNet
 - CMcsBus_AxisParametersNet, [198](#)
 - ~CMcsBus_ExtensionNet
 - CMcsBus_ExtensionNet, [199](#)
 - ~CMcsBus_FYIExtensionNet
 - CMcsBus_FYIExtensionNet, [200](#)
 - ~CMcsBus_MotorControlNet
 - CMcsBus_MotorControlNet, [205](#)
 - ~CMcsBus_SensorNet
 - CMcsBus_SensorNet, [221](#)
 - ~CMcsBus_TempSensorNet
 - CMcsBus_TempSensorNet, [230](#)
 - ~CMcsBus_VoltageModeNet
 - CMcsBus_VoltageModeNet, [232](#)
 - ~CMcsUsbDacqNet
 - CMcsUsbDacqNet, [246](#)
 - ~CMcsUsbFactoryNet
 - CMcsUsbFactoryNet, [291](#)
 - ~CMcsUsbFunctionNet
 - CMcsUsbFunctionNet, [299](#)
 - ~CMcsUsbListEntryNet
 - CMcsUsbListEntryNet, [300](#)
 - ~CMcsUsbListNet
 - CMcsUsbListNet, [306](#)
 - ~CMcsUsbNet
 - CMcsUsbNet, [312](#)
 - ~CMeFunctionNet
 - CMeFunctionNet, [363](#)
 - ~CMeaCleanDeviceNet
 - CMeaCleanDeviceNet, [335](#)
 - ~CMeaCoatDeviceNet
 - CMeaCoatDeviceNet, [339](#)
 - ~CMeaDeviceNet
 - CMeaDeviceNet, [345](#)
 - ~CMealImpedanceDeviceNet
 - CMealImpedanceDeviceNet, [358](#)
 - ~CMeaSwitchDeviceNet
 - CMeaSwitchDeviceNet, [360](#)
 - ~CMeaUSBDeviceNet
 - CMeaUSBDeviceNet, [362](#)
 - ~CMultiBatteryChargerDeviceNet
 - CMultiBatteryChargerDeviceNet, [365](#)
 - ~CMultiwellCallbackFunctionNet
 - CMultiwellCallbackFunctionNet, [373](#)
 - ~CMultiwellDeviceNet
 - CMultiwellDeviceNet, [376](#)
 - ~CMultiwellOptoStimFunctionNet
 - CMultiwellOptoStimFunctionNet, [383](#)
 - ~CNF_GenDeviceNet
 - CNF_GenDeviceNet, [387](#)
 - ~CokuvisionStimulatorDeviceNet
 - CokuvisionStimulatorDeviceNet, [392](#)
 - ~CPPCFunctionNet
 - CPPCFunctionNet, [416](#)
 - ~CPathIdentDeviceNet
 - CPathIdentDeviceNet, [396](#)
 - ~CPedoterDeviceNet
 - CPedoterDeviceNet, [397](#)
 - ~CPeristalticPumpDeviceNet
 - CPeristalticPumpDeviceNet, [399](#)
 - ~CPgaDeviceNet
 - CPgaDeviceNet, [400](#)
 - ~CPositionIIDeviceNet
 - CPositionIIDeviceNet, [403](#)
 - ~CPositionImpDeviceNet
 - CPositionImpDeviceNet, [411](#)
 - ~CProgramPressureCurveNet
 - CProgramPressureCurveNet, [429](#)
 - ~CPulseGeneratorFunctionNet
 - CPulseGeneratorFunctionNet, [431](#)
 - ~CRFFFunctionNet
 - CRFFFunctionNet, [439](#)
 - ~CRetinaLedDeviceNet
 - CRetinaLedDeviceNet, [437](#)
 - ~CRoboDeviceNet
 - CRoboDeviceNet, [464](#)
 - ~CRoboFluidDeviceNet
 - CRoboFluidDeviceNet, [478](#)
 - ~CSCUFunctionNet
 - CSCUFunctionNet, [495](#)
 - ~CSafeISDeviceNet
 - CSafeISDeviceNet, [489](#)
 - ~CStg200xBasicNet
 - CStg200xBasicNet, [513](#)
 - ~CStg200xDownloadNet
 - CStg200xDownloadNet, [551](#)
 - ~CSw2to64DeviceNet
 - CSw2to64DeviceNet, [568](#)
 - ~CTEERFunctionNet
 - CTEERFunctionNet, [585](#)
 - ~CTEERMachineDeviceNet
 - CTEERMachineDeviceNet, [593](#)
 - ~CTcxDeviceNet
 - CTcxDeviceNet, [572](#)
 - ~CUsbDeviceConfigurationFunctionNet
 - CUsbDeviceConfigurationFunctionNet, [595](#)
 - ~CWarnerUssingDeviceNet
 - CWarnerUssingDeviceNet, [612](#)
 - ~CWarnerUssingFunctionNet
 - CWarnerUssingFunctionNet, [615](#)
 - ~CWarnerValveControllerDeviceNet
 - CWarnerValveControllerDeviceNet, [634](#)
 - ~CWarnerValveControllerDeviceTesterFunctionNet
 - CWarnerValveControllerDeviceTesterFunctionNet, [655](#)
 - ~DriverVersionNet
 - DriverVersionNet, [668](#)
 - ~SidebandData
 - CStimulusFunctionNet::SidebandData, [694](#)
 - ~StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- A
- CFilterCoefficientsNet, [131](#)
 - Mcs::Usb, [67](#)
- AccelOnly

- Mcs::Usb, [87](#)
- AdapterTypeEnumNet
 - Mcs::Usb, [50](#)
- AdditionalText
 - HeadstageIDTypeObject, [681](#)
- AddLoopEntry
 - CRetinaLedDeviceNet, [437](#)
- AddSelectedChannelsQueue
 - CMcsUsbDacqNet, [246](#), [247](#), [249](#)
- AddSoftwareKey
 - CMcsUsbNet, [312](#)
- AddTableEntry
 - CRetinaLedDeviceNet, [437](#)
- ALA_VC3
 - Mcs::Usb, [73](#), [87](#)
- ALA_VC3_DEVICE
 - Mcs::Usb, [59](#)
- ALTERA
 - Mcs::Usb, [52](#)
- Altera
 - FirmwareDestinationNames, [674](#)
- ALTERA_BASE
 - Mcs::Usb, [53](#)
- ALTERA_BOOTSTRAP
 - Mcs::Usb, [53](#)
- ALTERA_GOLD
 - Mcs::Usb, [53](#)
- ALTERA_TARGET1
 - Mcs::Usb, [54](#)
- ALTERA_TARGET2
 - Mcs::Usb, [54](#)
- ALTERA_TARGET3
 - Mcs::Usb, [54](#)
- AlwaysOn
 - Mcs::Usb, [76](#)
- Ampere
 - Mcs::Usb, [51](#), [87](#)
- amplification
 - CMeaAudioFunctionNet::s_setaudionet, [693](#)
 - CW2100_FunctionNet::AudioChannelsNet, [92](#)
- AmplifierSettle
 - CIntanMea_FunctionNet, [183](#)
- AMS_Dongle
 - Mcs::Usb, [74](#)
- Analog
 - Mcs::Usb, [73](#), [90](#)
- AnalogGain
 - CMeaDeviceNet, [350](#)
- AnalogGroup
 - Mcs::Usb, [58](#)
- AnalogOut_DAC_Range_EnumNet
 - Mcs::Usb, [51](#)
- AnalogSource_HS1
 - Mcs::Usb, [51](#)
- AnalogSource_HS2
 - Mcs::Usb, [51](#)
- AnalogSource_IF
 - Mcs::Usb, [51](#)
- AnalogSourceEnumNet
 - Mcs::Usb, [51](#)
- AnalogUnitEnumNet
 - Mcs::Usb, [51](#)
- Any
 - Mcs::Usb, [73](#), [87](#)
- ApplyGains
 - CPgaDeviceNet, [400](#)
- AreTransistorVoltagesSet
 - CCMOSMea_FunctionNet, [99](#)
- Armed
 - Mcs::Usb, [81](#)
- ASMedia
 - Mcs::Usb, [86](#)
- AssociateToThis
 - CMcsUsbNet, [312](#)
- AudioTestChannelGroup
 - Mcs::Usb, [57](#), [69](#), [79](#), [89](#)
- AutomaticAnalogOut
 - CSCUFunctionNet, [495](#)
- Aux
 - Mcs::Usb, [60](#)
- AuxIn
 - Mcs::Usb, [60](#), [70](#), [79](#), [82](#), [89](#)
- AuxPort
 - Mcs::Usb, [56](#)
- Axes_I
 - CRoboDeviceNet, [472](#)
- Axes_X
 - CRoboDeviceNet, [472](#)
- Axes_Y
 - CRoboDeviceNet, [472](#)
- Axes_Z
 - CRoboDeviceNet, [473](#)
- Axis_I
 - CRoboDeviceNet, [473](#)
- Axis_X
 - CRoboDeviceNet, [473](#)
- Axis_Y
 - CRoboDeviceNet, [473](#)
- Axis_Z
 - CRoboDeviceNet, [473](#)
- B
 - CFilterCoefficientsNet, [131](#)
 - Mcs::Usb, [67](#)
- BatteryState, [92](#)
 - Charge, [92](#)
 - ChargeRegionString, [92](#)
 - ChargeString, [92](#)
 - Voltage, [93](#)
 - VoltageString, [93](#)
- BcdDevice
 - DeviceldNet, [664](#)
- BeginImpedanceCheck
 - CIntanMea_FunctionNet, [183](#)
- Bessel
 - Mcs::Usb, [64](#)
- BesselFilterHighPassNet, [93](#)

- BesselFilterHighPassNet, 93
- BesselFilterLowPassNet, 93
 - BesselFilterLowPassNet, 94
- BMI
 - Mcs::Usb, 83
- bmRequestType
 - usbSetupPacket_t, 697
- BOOST_BIT
 - CW2100_StimulatorFunctionNet, 611
- Bootstrap
 - FirmwareDestinationNames, 674
 - Mcs::Usb, 52
- BootstrapOtherCypress
 - Mcs::Usb, 52
- Both
 - Mcs::Usb, 87
- Break
 - Mcs::Usb, 77
- bRequest
 - usbSetupPacket_t, 697
- BurnAdcOffset
 - COctoPotDeviceNet, 388
- BurnDacOffset
 - CDacCalibrationFunctionNet, 113
 - COctoPotDeviceNet, 389
- BUS1_MCSBUS1
 - FirmwareDestinationNames, 674
- BUS1_MCSBUS2
 - FirmwareDestinationNames, 674
- BUS1MCSBUS1
 - Mcs::Usb, 52
- BUS1MCSBUS2
 - Mcs::Usb, 52
- BUS2MCSBUS1
 - Mcs::Usb, 52
- BUS2MCSBUS2
 - Mcs::Usb, 52
- BUSNUMBER1
 - Mcs::Usb, 52
- BUSNUMBER2
 - Mcs::Usb, 52
- BusType
 - DeviceIdNet, 664
- Butterworth
 - Mcs::Usb, 64
- ButterworthFilterHighPassNet, 94
 - ButterworthFilterHighPassNet, 94
- ButterworthFilterLowPassNet, 95
 - ButterworthFilterLowPassNet, 95
- CalibrateThermocouple
 - CFluidControlDeviceNet, 138
 - CTcxDeviceNet, 572
- Campden_Ci4600EphysVideoDataIntegrator
 - Mcs::Usb, 73
- CancelInternalCalibration
 - CTEERFunctionNet, 585
- CancelPoolLoop
 - CRoboDeviceNet, 464
- CancelPoolLoopAndStopMovement
 - CRoboDeviceNet, 464
- CancelTableLoop
 - CRoboDacqNet, 449
- CancelTableLoopAndStopTable
 - CRoboDacqNet, 449
- CapacityTest
 - CMultiBatteryChargerDeviceNet, 365
- CatchAmp
 - Mcs::Usb, 72
- CatchAmpGetAdcMean
 - CMcsBus_SensorNet, 221
- CatchAmpGetAdcValue
 - CMcsBus_SensorNet, 221
- CatchAmpGetAdcValueH
 - CMcsBus_SensorNet, 221
- CatchAmpGetAdcValueL
 - CMcsBus_SensorNet, 221
- CatchAmpGetDacAmplitude
 - CMcsBus_SensorNet, 222
- CatchAmpGetDacEnable
 - CMcsBus_SensorNet, 222
- CatchAmpGetDacOffset
 - CMcsBus_SensorNet, 222
- CatchAmpGetPwmEnable
 - CMcsBus_SensorNet, 222
- CatchAmpSetDacAmplitude
 - CMcsBus_SensorNet, 222
- CatchAmpSetDacEnable
 - CMcsBus_SensorNet, 222
- CatchAmpSetDacOffset
 - CMcsBus_SensorNet, 222
- CatchAmpSetPwmEnable
 - CMcsBus_SensorNet, 222
- CChannelTestDeviceNet, 95
 - ~CChannelTestDeviceNet, 96
 - CChannelTestDeviceNet, 96
 - SetAmplitude, 96
 - SetAttenuation, 96
 - SetFrequency, 96
 - SetWaveform, 96
- CCMOSMea_FunctionNet, 96
 - AreTransistorVoltagesSet, 99
 - CCMOSMea_FunctionNet, 98, 99
 - ClearSTGOutput, 99
 - DetectChipType, 99
 - EnableChannelsInGroup, 99
 - GetADCInputOffset, 99
 - GetBath, 99
 - GetBathMode, 99
 - GetEnabledChannelsInGroup, 100
 - GetGate, 100
 - GetGNDI, 100
 - GetGroupADCBits, 100
 - GetGroupChannelBitmaskBySelect, 100
 - GetGroupChannelBitmaskHS1NCBathCurrent, 100, 101

- GetGroupChannelBitmaskHS1NCCol2Current, 101
- GetGroupChannelBitmaskHS1NChipTemp, 101
- GetGroupChannelBitmaskHS1Sidebands, 101
- GetGroupChannelBitmaskHS1TriggerStatus, 101, 102
- GetGroupChannelBitmaskIFDigChannels, 102
- GetGroupChannelBitmaskInterfaceADC, 102
- GetGroupChannelBitmaskPacketFrameContext, 102
- GetGroupChannelBitmaskSTG1DACSignal, 102, 103
- GetGroupDCOffset, 103
- GetGroupID, 103
- GetGroupNumberOfChannels, 103
- GetGroupResolutionPerDigit, 103
- GetGroupSampleSize, 104
- GetGroupType, 104
- GetGroupUnit, 104
- GetMaxNumOfColumns, 104
- GetNeurochipMemoryData, 104
- GetNeurochipMemorySize, 105
- GetNumberOfSupportedGroups, 105
- GetSourceBulk, 105
- GetSourceDrain, 105
- GetSourceGate, 105
- GetStimulusSites, 105
- GetVDD3I, 105
- GetVDDI, 105
- IsChipPowered, 105
- IsGateFloating, 105
- PowerChip, 106
- SetADCInputOffset, 106
- SetBath, 106
- SetBathMode, 106
- SetGate, 106
- SetGateFloating, 106
- SetGateToVOP, 106
- SetNeurochipMemoryData, 106
- SetSourceBulk, 106
- SetSourceDrain, 107
- SetSourceGate, 107
- SetStimulusSites, 107
- UpdateTransistorVoltages, 107
- VOPSTimerSetResetTimes, 107
- CCMOSMeaDeviceNet, 107
 - ~CCMOSMeaDeviceNet, 108
 - CCMOSMeaDeviceNet, 108
 - CMosMea, 110
 - GetAvailableBaseSamplerates, 108
 - GetBaseSamplerate, 109
 - GetChannelDataI16, 109
 - GetChannelDataI32, 109
 - GetChannelDataUI16, 109
 - GetChannelDataUI32, 109
 - GetCMOSDataDictionary, 109
 - GetMaxReadableColumns, 109
 - SetBaseSamplerate, 110
 - SetRegionOfInterests, 110
 - Stimulus, 110
 - UpdateChannelBlock, 110
- CCMOSMeaDeviceNet::CRegionOfInterestRect, 435
 - CRegionOfInterestRect, 435
 - DeepCopy, 435
 - m_Bottom, 436
 - m_Left, 436
 - m_Right, 436
 - m_Top, 436
- CCreateFilterNet, 110
 - ~CCreateFilterNet, 111
 - CCreateFilterNet, 111
 - CutoffFrequency, 112
 - FindFilter, 111
 - GetBiQuad, 112
 - GetBiQuads, 112
 - NumCoefSets, 112
 - Order, 112
 - SampleRate, 112
 - Scale, 112
- CDacCalibrationFunctionNet, 112
 - !CDacCalibrationFunctionNet, 113
 - ~CDacCalibrationFunctionNet, 113
 - BurnDacOffset, 113
 - CDacCalibrationFunctionNet, 113
 - GetDacOffset, 114
 - SetDacOffset, 114
- CDacqGroupChannelGenericSelectionNet, 114
 - CDacqGroupChannelGenericSelectionNet, 115
- CDacqGroupChannelSelectionNet, 115
 - CDacqGroupChannelSelectionNet, 115
- CDacqGroupChannelSelectionTemplateNet
 - CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >, 116
- CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >, 115
 - CDacqGroupChannelSelectionTemplateNet, 116
 - EnableChannelsInGroup, 116
 - GetDeviceGroupChannelInfos, 116, 117
 - GetEnabledChannelsInGroup, 117
 - GetGroupID, 117
 - GetGroupNumberOfChannels, 117
 - GetGroupSampleSize, 117
 - GetGroupType, 118
 - GetNumberOfSupportedGroups, 118
- CDeviceGroupChannelInfoGenericNet, 118
 - CDeviceGroupChannelInfoGenericNet, 118
- CDeviceGroupChannelInfoMEA2100_256Net, 119
 - CDeviceGroupChannelInfoMEA2100_256Net, 119
- CDeviceGroupChannelInfoNet, 119
 - CDeviceGroupChannelInfoNet, 120
- CDeviceGroupChannelInfoSCUNet, 120
 - CDeviceGroupChannelInfoSCUNet, 120

- CDeviceGroupChannelInfoTemplateNet
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, 121
- CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, 120
 - CDeviceGroupChannelInfoTemplateNet, 121
 - GroupID, 121
 - GroupType, 121
 - NumberOfChannels, 121
- CDeviceGroupChannelInfoW2100Net, 121
 - CDeviceGroupChannelInfoW2100Net, 122
- CDigOutStimulatorFunctionNet, 122
 - !CDigOutStimulatorFunctionNet, 123
 - ~CDigOutStimulatorFunctionNet, 123
 - CDigOutStimulatorFunctionNet, 123
 - ClearChannel, 123
 - GetGlobalRepeat, 123
 - GetNumberOfChannels, 124
 - GetStartTriggerSlope, 124
 - GetStopTriggerSlope, 124
 - PrepareChannelData, 125
 - SendChannelData, 125
 - SetGlobalRepeat, 125
 - SetStartTriggerSlope, 126
 - SetStopTriggerSlope, 126
- CEncapsulatorDeviceNet, 126
 - CEncapsulatorDeviceNet, 127
 - GetRoboFluidDevice, 127
- CExternDTesterDeviceNet, 127
 - !CExternDTesterDeviceNet, 128
 - ~CExternDTesterDeviceNet, 128
 - CExternDTesterDeviceNet, 128
 - Read, 128
 - Read2, 128
 - Write, 128
 - Write2, 129
- CFilterCoefficientsNet, 129
 - ~CFilterCoefficientsNet, 130
 - A, 131
 - B, 131
 - CFilterCoefficientsNet, 129, 130
 - GetUIntA, 130
 - GetUIntB, 130
 - IsEqual, 130
- CFilterCoefficientsNet::s_FilterAttributesNet, 692
 - CommaPositionA, 693
 - CommaPositionB, 693
 - PostCommaA, 693
 - PostCommaB, 693
 - PreCommaA, 693
 - PreCommaB, 693
 - s_FilterAttributesNet, 692
 - ToCpp, 692
- CFilterConfigurationNet, 131
 - CFilterConfigurationNet, 131
 - EraseFilterParameterPermanent, 132
 - GetFilterAttributes, 132
 - GetHighpassFilterEnable, 132
 - ResetHighpassFilter, 132
 - SetFilterParameter, 132
 - SetFilterParameterPermanent, 132
 - SetHighpassFilterEnable, 133
- CFilterConfigurationRegisterNet, 133
 - CFilterConfigurationRegisterNet, 133
 - EraseFilterParameterPermanent, 133, 134
 - SetFilterParameter, 134
 - SetFilterParameterPermanent, 134
- CFilterPropertyNet, 135
 - ~CFilterPropertyNet, 135
 - CFilterPropertyNet, 135
 - CornerFrequency, 135
 - CornerFrequencymHz, 136
 - FilterActive, 136
 - FilterBand, 136
 - FilterFamily, 136
 - FilterType, 136
 - Order, 136
 - ToString, 135
- CFirmwareDestinationNet
 - Mcs::Usb, 51
- CFluidControlDeviceNet, 136
 - ~CFluidControlDeviceNet, 138
 - CalibrateThermocouple, 138
 - CFluidControlDeviceNet, 138
 - GetAdc, 138
 - GetDigin, 138
 - GetDigout, 139
 - GetPWM, 139
 - GetReferenceTemperature, 139
 - GetSingleValve, 139
 - GetThermocoupleCalibration, 140
 - GetThermocoupleNanovoltPerKelvin, 140
 - GetThermocoupleTemperature, 140
 - GetValve, 141
 - McsBus_VoltageMode, 143
 - SetDigout, 141
 - SetPWM, 141
 - SetSingleValve, 141
 - SetThermocoupleNanovoltPerKelvin, 143
 - SetValve, 143
- CFYIDeviceNet, 143
 - CFYIDeviceNet, 144
 - FYIProgram, 144
 - FYITemp, 144
 - Sensor, 144
- CGenericDevelopDeviceNet, 145
 - ~CGenericDevelopDeviceNet, 152
 - CGenericDevelopDeviceNet, 152
 - ClosePipe, 152
 - FindEndpoints, 152
 - GetBuffer, 152
 - GetByteBuffer, 153
 - GetIntBuffer, 153
 - GetShortBuffer, 154
 - GetUByteBuffer, 155
 - GetUIntBuffer, 155

- GetUShortBuffer, 156
- OpenPipe, 157
- ReadPipe, 157
- ResetPipe, 157
- SetBuffer, 158
- SetByteBuffer, 158
- SetIntBuffer, 158
- SetShortBuffer, 159
- SetUByteBuffer, 160
- SetUIntBuffer, 160
- SetUShortBuffer, 161
- SetValue, 161
- VendorInRequest, 162
- VendorOutRequest, 162
- WritePipe, 162
- CGilsonDeviceNet, 163
 - ~CGilsonDeviceNet, 164
 - CGilsonDeviceNet, 163
 - ConnectSlave, 164
 - GetLastAnswer, 164
 - m_pGilsonDevice, 164
 - SendBuffered, 164
 - SendImmediate, 164
 - SendImmediateGetResponse, 164
- CGrapheneASICDeviceNet, 165
 - ~CGrapheneASICDeviceNet, 165
 - CGrapheneASICDeviceNet, 165
 - GetAvailableBaseSamplerates, 165
 - SetBaseSamplerate, 166
 - SetRegionOfInterests, 166
- CGrapheneFunctionNet, 166
 - !CGrapheneFunctionNet, 168
 - ~CGrapheneFunctionNet, 168
 - CGrapheneFunctionNet, 168
 - GetCur2VolResistance, 168
 - GetDACOffset, 169
 - GetVdsVgs, 169, 170
 - GetVdVs, 170
 - GetVdVsDAC, 170, 172
 - GetVoltageRange, 172
 - GetVoltageReached, 172, 173
 - GetVoltageResolution, 173
 - SetDACOffset, 173, 174
 - SetVdsVgs, 174
 - SetVdVs, 175
 - SetVdVsDAC, 175
 - SetVoltageRange, 176
 - SetVoltageResolution, 176
- ChangeSerialNumber
 - CMcsUsbFactoryNet, 291
- channel
 - CMeaAudioFunctionNet::s_setaudionet, 693
 - CW2100_FunctionNet::AudioChannelsNet, 92
- ChannelBlock_AvailFrames
 - CMcsUsbDacqNet, 249
- ChannelBlock_ReadAsFrameArrayI16
 - CMcsUsbDacqNet, 249, 250
- ChannelBlock_ReadAsFrameArrayI32
 - CMcsUsbDacqNet, 251
- ChannelBlock_ReadAsFrameArrayUI16
 - CMcsUsbDacqNet, 252, 253
- ChannelBlock_ReadAsFrameArrayUI32
 - CMcsUsbDacqNet, 253, 254
- ChannelBlock_ReadFramesDictI16
 - CMcsUsbDacqNet, 255
- ChannelBlock_ReadFramesDictI32
 - CMcsUsbDacqNet, 255
- ChannelBlock_ReadFramesDictUI16
 - CMcsUsbDacqNet, 256
- ChannelBlock_ReadFramesDictUI32
 - CMcsUsbDacqNet, 257
- ChannelBlock_ReadFramesI16
 - CMcsUsbDacqNet, 257, 258
- ChannelBlock_ReadFramesI32
 - CMcsUsbDacqNet, 258, 259
- ChannelBlock_ReadFramesUI16
 - CMcsUsbDacqNet, 260
- ChannelBlock_ReadFramesUI32
 - CMcsUsbDacqNet, 261, 262
- channeldata_current
 - Mcs::Usb, 82
- channeldata_current_always_boost
 - Mcs::Usb, 82
- channeldata_current_always_boost_own_sync
 - Mcs::Usb, 82
- channeldata_current_own_boost_gnd_sync
 - Mcs::Usb, 82
- channeldata_current_own_sync
 - Mcs::Usb, 82
- channeldata_positive_current
 - Mcs::Usb, 82
- channeldata_positive_current_own_boost_gnd_sync
 - Mcs::Usb, 82
- channeldata_positive_current_own_sync
 - Mcs::Usb, 82
- channeldata_positive_voltage
 - Mcs::Usb, 82
- channeldata_voltage
 - Mcs::Usb, 82
- ChannelDataEvent
 - CMcsUsbDacqNet, 287
- ChannelPIC
 - Mcs::Usb, 52
- ChannelReset
 - CMultiBatteryChargerDeviceNet, 366
- ChannelTest
 - Mcs::Usb, 75
- Charge
 - BatteryState, 92
- ChargeRegionString
 - BatteryState, 92
- ChargeString
 - BatteryState, 92
- ChecksumAndPacketCounter
 - Mcs::Usb, 56
- CHiClampDeviceNet, 177

- CHiClampDeviceNet, 177
- RoboDacq, 178
- CHLADacqNet, 178
 - CHLADacqNet, 178
- CHLADeviceNet, 178
 - CHLADeviceNet, 179
 - HLADacq, 179
 - SerialPort, 179
- CHWInfo
 - CMcsUsbDacqNet::CHWInfo, 180
- Ci4600Intan
 - Mcs::Usb, 50
- CIntanMea_FunctionNet, 182
 - AmplifierSettle, 183
 - BeginImpedanceCheck, 183
 - CIntanMea_FunctionNet, 183
 - GetDSPHighPassByIndex, 183
 - GetImpedanceResult, 183
 - GetIntanRegister, 183
 - GetLowerFrequencyByIndex, 183
 - GetUpperFrequencyByIndex, 184
 - SetBandwidthByIndex, 184
 - SetDiagnosticMode, 184
 - SetDSPHighPassByIndex, 184
 - SetIntanRegister, 184
- CInterfaceboard2FunctionNet, 184
 - !CInterfaceboard2FunctionNet, 185
 - ~CInterfaceboard2FunctionNet, 185
 - CInterfaceboard2FunctionNet, 185
 - GetIoVoltage, 185
 - SetIoVoltage, 186
- CInterfaceboardFunctionNet, 186
 - !CInterfaceboardFunctionNet, 187
 - ~CInterfaceboardFunctionNet, 187
 - CInterfaceboardFunctionNet, 187
 - GetCardinalDacqSamplerate, 187
 - GetCardinalStgOutputrate, 187
 - SetCardinalDacqSamplerate, 187
 - SetCardinalStgOutputrate, 188
- ClampAmpRestart
 - CRoboDacqNet, 449
- ClampModeCurrent
 - Mcs::Usb, 84
- ClampModeInternalCalibration
 - Mcs::Usb, 84
- ClampModeOpen
 - Mcs::Usb, 84
- ClampModeVoltage
 - Mcs::Usb, 83
- ClearBuffers
 - CMcsUsbDacqNet, 262
- ClearChannel
 - CDigOutStimulatorFunctionNet, 123
- ClearChannel_PrepAndSendData
 - CStg200xDownloadNet, 551
 - CStimulusFunctionNet, 558
- ClearChannelData
 - CStg200xDownloadBasicNet, 543
- CStimulusFunctionNet, 558
 - CW2100_StimulatorFunctionNet, 606
- ClearMultiplexedData
 - CStimulusFunctionNet, 558
- ClearSTGOutput
 - CCMOSMea_FunctionNet, 99
- ClearStimulusParametersCache
 - CW2100_FunctionNet, 599
- ClearSyncData
 - CStg200xDownloadBasicNet, 544
 - CStimulusFunctionNet, 558
- ClearTable
 - CRetinaLedDeviceNet, 437
- ClearTableName
 - CWarnerValveControllerDeviceNet, 634
- ClearUserDefinedNameCache
 - CW2100_FunctionNet, 599
- ClearValveTable
 - CWarnerValveControllerDeviceNet, 634
- CLIH3DeviceNet, 188
 - !CLIH3DeviceNet, 190
 - ~CLIH3DeviceNet, 190
 - CLIH3DeviceNet, 190
 - DummyCommand, 190
 - EnableUserTrigger, 190
 - ErasePermanentAdcOffset, 191
 - ErasePermanentDacOffset, 191
 - GetAdcOffset, 191
 - GetAudioOutDacParameter, 191
 - GetDacIdleValue, 192
 - GetDacOffset, 192
 - GetDacqRunStatus, 192
 - GetDacUseldleValue, 192
 - GetDigInState, 193
 - GetEEPROMPage, 193
 - GetSampleInterval, 193
 - IsUserTriggerEnabled, 193
 - ReadClipping, 194
 - ReadUARTData, 194
 - SendCommand, 194
 - SetAdcOffset, 194
 - SetAdcOffsetPermanent, 195
 - SetAudioOutDacParameter, 195
 - SetDacIdleValue, 195
 - SetDacOffset, 195
 - SetDacOffsetPermanent, 196
 - SetDacUseldleValue, 196
 - SetDigOutState, 196
 - SetEEPROMPage, 196
 - SetSampleInterval, 197
 - StimulusFunction, 197
 - WriteUARTData, 197
- Close
 - Mcs::Usb, 72
- CloseAllValves
 - CRoboFluidDeviceNet, 478
- ClosePipe
 - CGenericDevelopDeviceNet, 152

ClosePlateClamp
 CMultiwellDeviceNet, 376

CMcsBus_AxisParametersNet, 197
 ~CMcsBus_AxisParametersNet, 198
 CMcsBus_AxisParametersNet, 198
 GetAxisParametersSignedEeprom, 198
 GetAxisParametersUnsignedEeprom, 198
 SetAxisParametersEeprom, 198, 199

CMcsBus_ExtensionNet, 199
 ~CMcsBus_ExtensionNet, 199
 CMcsBus_ExtensionNet, 199
 GetLEDSwitch, 200
 SetLEDSwitch, 200

CMcsBus_FYIExtensionNet, 200
 ~CMcsBus_FYIExtensionNet, 200
 CMcsBus_FYIExtensionNet, 200
 GetDIO, 201
 GetSingleHeater, 201
 GetValves, 201
 SetDIO, 201
 SetSingleHeater, 201
 SetValves, 201

CMcsBus_MotorControlNet, 202
 ~CMcsBus_MotorControlNet, 205
 CMcsBus_MotorControlNet, 205
 GetMCAcceleration, 205
 GetMCAccelerationEeprom, 205
 GetMCAccelerationShortCommand, 205
 GetMCAXisRevisionEeprom, 205
 GetMCBreakCurrent, 205
 GetMCBreakCurrentEeprom, 206
 GetMCConfig, 206
 GetMCConfigEeprom, 206
 GetMCCurrent, 206
 GetMCCurrentEeprom, 206
 GetMCCurrentMode, 206
 GetMCCurrentModeEeprom, 206
 GetMCCurrentModeShortCommand, 207
 GetMCCurrentPosition, 207
 GetMCCurrentShortCommand, 207
 GetMCCurrentSpeed, 207
 GetMCMaxAcceleration, 207
 GetMCMaxAccelerationEeprom, 207
 GetMCMaxCurrent, 207
 GetMCMaxCurrentEeprom, 208
 GetMCMaxSpeed, 208
 GetMCMaxSpeedEeprom, 208
 GetMCMaxTravel, 208
 GetMCMaxTravelEeprom, 208
 GetMCMaxTravelShortCommand, 208
 GetMCMovement, 208
 GetMCNewPosition, 209
 GetMCOOutputOnOff, 209
 GetMCPhase, 209
 GetMCPhaseOffset, 209
 GetMCReference, 209
 GetMCReferenceCurrent, 209
 GetMCReferenceCurrentEeprom, 209
 GetMCRegulatorGain, 210
 GetMCRegulatorGainEeprom, 210
 GetMCScalingFactor, 210
 GetMCScalingFactorEeprom, 210
 GetMCSpeed, 210
 GetMCSpeedEeprom, 210
 GetMCSpeedShortCommand, 210
 GetMCSpeedUnitEeprom, 211
 GetMCStandbyCurrent, 211
 GetMCStandbyCurrentEeprom, 211
 GetMCStandbyTime, 211
 GetMCStandbyTimeEeprom, 211
 GetSubChannel, 211
 SetMCAcceleration, 211
 SetMCAccelerationEeprom, 212
 SetMCAccelerationShortCommand, 212
 SetMCAXisRevisionEeprom, 212
 SetMCBreakCurrent, 212
 SetMCBreakCurrentEeprom, 212
 SetMCConfig, 212
 SetMCConfigEeprom, 213
 SetMCCurrent, 213
 SetMCCurrentEeprom, 213
 SetMCCurrentMode, 213
 SetMCCurrentModeEeprom, 213
 SetMCCurrentModeShortCommand, 213
 SetMCCurrentPosition, 214
 SetMCCurrentShortCommand, 214
 SetMCMaxAcceleration, 214
 SetMCMaxAccelerationEeprom, 214
 SetMCMaxCurrent, 214
 SetMCMaxCurrentEeprom, 214
 SetMCMaxSpeed, 215
 SetMCMaxSpeedEeprom, 215
 SetMCMaxTravel, 215
 SetMCMaxTravelEeprom, 215
 SetMCMaxTravelShortCommand, 215
 SetMCNewPosition, 215
 SetMCOOutputOnOff, 216
 SetMCReference, 216
 SetMCReferenceCurrent, 216
 SetMCReferenceCurrentEeprom, 216
 SetMCRegulatorGain, 216
 SetMCRegulatorGainEeprom, 216
 SetMCRotation, 217
 SetMCScalingFactor, 217
 SetMCScalingFactorEeprom, 217
 SetMCSpeed, 217
 SetMCSpeedEeprom, 217
 SetMCSpeedShortCommand, 217
 SetMCSpeedUnitEeprom, 218
 SetMCStandbyCurrent, 218
 SetMCStandbyCurrentEeprom, 218
 SetMCStandbyTime, 218
 SetMCStandbyTimeEeprom, 218
 SetSubChannel, 218
 StartMCMovement, 219
 StopMCMovement, 219

- CMcsBus_SensorNet, [219](#)
 - ~CMcsBus_SensorNet, [221](#)
 - CatchAmpGetAdcMean, [221](#)
 - CatchAmpGetAdcValue, [221](#)
 - CatchAmpGetAdcValueH, [221](#)
 - CatchAmpGetAdcValueL, [221](#)
 - CatchAmpGetDacAmplitude, [222](#)
 - CatchAmpGetDacEnable, [222](#)
 - CatchAmpGetDacOffset, [222](#)
 - CatchAmpGetPwmEnable, [222](#)
 - CatchAmpSetDacAmplitude, [222](#)
 - CatchAmpSetDacEnable, [222](#)
 - CatchAmpSetDacOffset, [222](#)
 - CatchAmpSetPwmEnable, [222](#)
 - CMcsBus_SensorNet, [221](#)
 - Get2AnalogInput, [223](#)
 - Get2DigitalInput, [223](#)
 - Get4ADC, [223](#)
 - Get4ADCAverage, [223](#)
 - Get4ADCCatchampAverageShift, [223](#)
 - Get4ADCMode, [223](#)
 - Get4DAC, [223](#)
 - GetADCs, [223](#)
 - GetADCsLoop, [224](#)
 - GetBubbleStatus, [224](#)
 - GetDACs, [224](#)
 - GetDetectionThreshold, [224](#)
 - GetDetectorValue, [224](#)
 - GetLatency, [224](#)
 - GetLatencyCounter, [224](#)
 - GetMinimalThreshold, [224](#)
 - GetMovePump, [225](#)
 - GetPiezoState, [225](#)
 - GetPressure, [225](#)
 - GetPressureOffset, [225](#)
 - GetRegulationTimeouts, [225](#)
 - GetRegulatorFactor, [226](#)
 - GetRegulatorOnOff, [226](#)
 - GetRegulatorStatus, [226](#)
 - GetRotatePump, [226](#)
 - GetSamplePeriode, [226](#)
 - GetSollPressure, [226](#)
 - GetSyncState, [226](#)
 - Set4ADCCatchampAverageShift, [227](#)
 - Set4ADCMode, [227](#)
 - Set4DAC, [227](#)
 - SetDACs, [227](#)
 - SetDetectionThreshold, [227](#)
 - SetLatency, [227](#)
 - SetMinimalThreshold, [227](#)
 - SetMovePump, [228](#)
 - SetPiezoState, [228](#)
 - SetPressureOffset, [228](#)
 - SetRegulationTimeouts, [228](#)
 - SetRegulatorFactor, [228](#)
 - SetRegulatorOnOff, [228](#)
 - SetRotatePump, [228](#)
 - SetSamplePeriode, [229](#)
 - SetSollPressure, [229](#)
 - StartSync, [229](#)
 - TactSwitchGetState, [229](#)
 - TactSwitchSetDisplay, [229](#)
- CMcsBus_TempSensorNet, [229](#)
 - ~CMcsBus_TempSensorNet, [230](#)
 - CMcsBus_TempSensorNet, [230](#)
 - GetNanoVoltsPerKelvin, [230](#)
 - GetTemperatur, [230](#)
 - GetThermoOffset, [230](#)
 - GetThermoTemp, [231](#)
 - GetThermoVoltage, [231](#)
 - SetNanoVoltsPerKelvin, [231](#)
 - SetThermoOffset, [231](#)
- CMcsBus_VoltageModeNet, [231](#)
 - ~CMcsBus_VoltageModeNet, [232](#)
 - CMcsBus_VoltageModeNet, [232](#)
 - GetVMMaxNegativeCurrent, [233](#)
 - GetVMMaxNegativeCurrentEeprom, [233](#)
 - GetVMMaxNegativeVoltage, [233](#)
 - GetVMMaxNegativeVoltageEeprom, [233](#)
 - GetVMMaxPositiveCurrent, [233](#)
 - GetVMMaxPositiveCurrentEeprom, [233](#)
 - GetVMMaxPositiveVoltage, [233](#)
 - GetVMMaxPositiveVoltageEeprom, [234](#)
 - GetVMOutputOnOff, [234](#)
 - GetVMVoltage, [234](#)
 - SetVMMaxNegativeCurrent, [234](#)
 - SetVMMaxNegativeCurrentEeprom, [234](#)
 - SetVMMaxNegativeVoltage, [234](#)
 - SetVMMaxNegativeVoltageEeprom, [234](#)
 - SetVMMaxPositiveCurrent, [235](#)
 - SetVMMaxPositiveCurrentEeprom, [235](#)
 - SetVMMaxPositiveVoltage, [235](#)
 - SetVMMaxPositiveVoltageEeprom, [235](#)
 - SetVMOutputOnOff, [235](#)
 - SetVMVoltage, [235](#)
- CMcsBusNet, [236](#)
 - ~CMcsBusNet, [236](#)
 - CMcsBusNet, [236](#)
 - CMcsBusNet::GetMode, [237](#)
 - CMcsBusNet::GetModeEeprom, [237](#)
 - CMcsBusNet::SetMode, [237](#)
 - CMcsBusNet::SetModeEeprom, [237](#)
 - GetBusAddress, [237](#)
 - GetBusAddressEeprom, [237](#)
 - GetCommand, [237, 238](#)
 - GetHWRevisionEeprom, [238](#)
 - SetBusAddress, [238](#)
 - SetBusAddressEeprom, [238](#)
 - SetCommand, [239](#)
 - SetHWRevisionEeprom, [239](#)
- CMcsBusNet::GetMode
 - CMcsBusNet, [237](#)
- CMcsBusNet::GetModeEeprom
 - CMcsBusNet, [237](#)
- CMcsBusNet::SetMode
 - CMcsBusNet, [237](#)

- CMcsBusNet::SetModeEeprom
 - CMcsBusNet, [237](#)
- CMcsUsbDacqNet, [240](#)
 - ~CMcsUsbDacqNet, [246](#)
 - AddSelectedChannelsQueue, [246](#), [247](#), [249](#)
 - ChannelBlock_AvailFrames, [249](#)
 - ChannelBlock_ReadAsFrameArrayI16, [249](#), [250](#)
 - ChannelBlock_ReadAsFrameArrayI32, [251](#)
 - ChannelBlock_ReadAsFrameArrayUI16, [252](#), [253](#)
 - ChannelBlock_ReadAsFrameArrayUI32, [253](#), [254](#)
 - ChannelBlock_ReadFramesDictI16, [255](#)
 - ChannelBlock_ReadFramesDictI32, [255](#)
 - ChannelBlock_ReadFramesDictUI16, [256](#)
 - ChannelBlock_ReadFramesDictUI32, [257](#)
 - ChannelBlock_ReadFramesI16, [257](#), [258](#)
 - ChannelBlock_ReadFramesI32, [258](#), [259](#)
 - ChannelBlock_ReadFramesUI16, [260](#)
 - ChannelBlock_ReadFramesUI32, [261](#), [262](#)
 - ChannelDataEvent, [287](#)
 - ClearBuffers, [262](#)
 - CMcsUsbDacqNet, [246](#)
 - CMcsUsbDacqNet::GetFilterProperties, [262](#)
 - Error_Callback_Aquisition_Stopped, [286](#)
 - Error_Callback_Data_Lost, [286](#)
 - Error_Callback_Frames_Lost, [286](#)
 - Error_Callback_Packet_Error, [286](#)
 - Error_Callback_Queue_Full, [287](#)
 - Error_Callback_RingQueue_Full, [287](#)
 - ErrorEvent, [287](#)
 - GetAdapterType, [262](#)
 - GetAdcDataFormat, [262](#)
 - GetAdcZero, [263](#)
 - GetAnalogValueUnit, [263](#)
 - GetChannelDataFillSize, [263](#)
 - GetChannelLayout, [263](#)
 - GetChannelsInBlock, [263](#)
 - GetDataFormat, [263](#)
 - GetDataMode, [263](#)
 - GetDigitalSource, [264](#), [265](#)
 - GetErrorMessage, [266](#)
 - GetFilterProperty, [266](#)
 - GetGroupChannelDataI16, [266](#)
 - GetGroupChannelDataI32, [266](#)
 - GetGroupChannelDataUI16, [267](#)
 - GetGroupChannelDataUI32, [267](#)
 - GetHardwareMaxRange, [268](#)
 - GetHardwareMinRange, [268](#)
 - GetMaxSamplingFrequency, [268](#)
 - GetMeaLayout, [268](#)
 - GetMinSamplingFrequencyStepsize, [269](#)
 - GetNumberOfDataBits, [269](#)
 - GetPoti, [269](#)
 - GetResolutionPerDigit, [269](#)
 - GetSamplerate, [269](#)
 - GetVoltageRangeIndex, [269](#)
 - GetVoltageRangeInMicroVolt, [270](#)
 - GetVoltageRangeInMilliVolt, [270](#)
 - HWInfo, [270](#)
 - Samplerate, [287](#)
 - SendStartDacq, [270](#)
 - SendStartStgAndDacq, [270](#)
 - SendStopDacq, [271](#)
 - SendStopStgAndDacq, [271](#)
 - SendStopStgAndDacqWithOptions, [271](#)
 - SetDataMode, [272](#)
 - SetDigitalSource, [272–274](#)
 - SetPoti, [274](#)
 - SetSamplerate, [274](#)
 - SetSelectedChannels, [275–277](#)
 - SetSelectedChannelsQueue, [277–279](#)
 - SetSelectedData, [279–281](#)
 - SetupGroupDacqQueue, [281](#)
 - SetVoltageRangeByIndex, [281](#)
 - SetVoltageRangeInMicroVolt, [282](#)
 - StartDacq, [282](#), [283](#)
 - StartLoop, [284](#), [285](#)
 - StopDacq, [286](#)
 - StopLoop, [286](#)
- CMcsUsbDacqNet::CHWInfo, [179](#)
 - CHWInfo, [180](#)
 - GetAvailableSampleRates, [180](#)
 - GetAvailableVoltageRangesInMicroVolt, [180](#)
 - GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt, [181](#)
 - GetNumberOfHWADCCChannels, [181](#)
 - GetNumberOfHWDigitalChannels, [181](#)
 - IsDigitalChannelDedicated, [182](#)
- CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, [596](#)
 - CVoltageRangeInfoNet, [597](#)
 - VoltageRangeDisplayStringMilliVolt, [597](#)
 - VoltageRangeInMicroVolt, [597](#)
- CMcsUsbDacqNet::GetFilterProperties
 - CMcsUsbDacqNet, [262](#)
- CMcsUsbDeviceStatePushFunctionNet, [287](#)
 - CMcsUsbDeviceStatePushFunctionNet, [288](#)
 - McsUsbDeviceStateEvent, [288](#)
 - TriggerStatus, [288](#)
- CMcsUsbDeviceStatePushNet, [288](#)
 - CMcsUsbDeviceStatePushNet, [289](#)
 - McsUsbDeviceStateEvent, [289](#)
 - TriggerStatus, [289](#)
- CMcsUsbFactoryNet, [289](#)
 - ~CMcsUsbFactoryNet, [291](#)
 - ChangeSerialNumber, [291](#)
 - CMcsUsbFactoryNet, [291](#)
 - Coldstart, [291](#)
 - CompareFirmware, [292](#)
 - DownloadFirmware, [292](#)
 - FindFirmwareVersionMagicInBuffer, [292](#)
 - FX3MCSDDataAddress, [297](#)
 - FX3MCSDDataDeviceIdOffset, [297](#)
 - FX3MCSDDataIFB1ImageOffset, [297](#)
 - FX3MCSDDataIFB2ImageOffset, [297](#)
 - FX3MCSDDataVersionOffset, [297](#)
 - GetChecksumFromFX3Image, [292](#)

- GetDestination, 292
- GetDestinationDisplayLabel, 292
- GetDestinationName, 292, 293
- GetDestinationSerialNumber, 293
- GetDestinationTargetAddress, 293
- GetFirmwareVersionFromFile, 293
- GetFirmwareVersionFromHexFile, 293
- GetNumDestinations, 293
- GetUSBDeviceIDFromFX3Image, 294
- GetUserCodeFromBitFile, 294
- GetUserCodeFromFlash, 294
- GetXilinxFlashOffset, 294
- GetXilinxFlashReadCommand, 294
- LoadUserFirmware, 294, 295
- ReadBlockFromFlash, 295
- ReadBlockFromIFBGlobalEEPROM, 295
- ReadBlockFromNVMEM, 295
- SetDestinationSerialNumber, 295
- UpdateFirmware, 295–297
- CMcsUsbFunctionNet, 298
 - !CMcsUsbFunctionNet, 299
 - ~CMcsUsbFunctionNet, 299
 - CMcsUsbFunctionNet, 298, 299
 - m_pMcsUsb, 299
 - m_pMcsUsbFunction, 299
 - ThrowCUsbExceptionNetOnError, 299
- CMcsUsbFunctionPointerContainer, 299
- CMcsUsbListEntryNet, 299
 - ~CMcsUsbListEntryNet, 300
 - DeviceId, 304
 - DeviceName, 304
 - Equals, 301
 - GetEntry, 302
 - GetEntryCount, 303
 - HwVersion, 304
 - Manufacturer, 304
 - Product, 304
 - SerialNumber, 304
 - SetStringFormat, 303
 - ToString, 304
- CMcsUsbListNet, 305
 - !CMcsUsbListNet, 306
 - ~CMcsUsbListNet, 306
 - CMcsUsbListNet, 305
 - Count, 307
 - DeviceArrival, 307
 - DeviceRemoval, 307
 - GetNumberOfDevices, 306
 - GetUsbListEntries, 306
 - GetUsbListEntry, 306
 - IsDeviceTypeOf, 307
 - SetStringFormat, 307
- CMcsUsbNet, 308
 - !CMcsUsbNet, 312
 - ~CMcsUsbNet, 312
 - AddSoftwareKey, 312
 - AssociateToThis, 312
 - CMcsUsbNet, 312
 - Connect, 312, 313
 - CyclePort, 314
 - Disconnect, 314
 - EmptyKey, 314
 - EnableExceptions, 314
 - EraseEEPROMRegisterPreconfig, 314
 - GetConfiguration, 315
 - GetDeviceCannotStallOutRequests, 315
 - GetDeviceCapableSpeed, 315
 - GetDeviceEnum, 315
 - GetDeviceId, 315
 - GetDeviceRootHubVendorEnum, 315
 - GetDeviceRootHubVendorID, 315
 - GetDeviceRootHubVendorName, 315
 - GetDeviceSpeed, 316
 - GetErrorText, 316
 - GetFirmwareVersion, 316
 - GetHardwareRevision, 316
 - GetHeadstageActive, 317
 - GetHeadstageID, 317
 - GetHeadstagePresent, 317
 - GetLastUSBError, 317
 - GetMea21UsbPort, 318
 - GetNumConfigurations, 318
 - GetSerialNumber, 318
 - GetSoftwareKey, 318
 - GetSoftwareKeyString, 318
 - GetStatus, 318
 - GetStatusOfLastCommand, 318
 - GetUsbListEntry, 319
 - GetVersion, 319
 - HasSoftwareKey, 319
 - IsConnected, 319
 - IsDeviceHighSpeed, 319
 - IsDeviceHighSpeedCapable, 319
 - IsExceptionsEnabled, 320
 - MultibootGetCypressImageId, 320
 - MultibootGetImageId, 320
 - MultibootGetSelectedImage, 320
 - MultibootSelectImage, 320
 - ReadEEPROMRegisterPreconfig, 320, 321
 - ReadRegister, 321
 - ReadRegister32, 321
 - ReadRegisterTimeSlot, 321
 - RemoveSoftwareKey, 321
 - RescanHeadstage, 321
 - SerialNumber, 328
 - SetConfiguration, 322
 - SetSoftwareKey, 322
 - Status_AlreadyConfigured, 324
 - Status_BadStartFrame, 324
 - Status_Btstuff, 324
 - Status_BufferOverrun, 324
 - Status_BufferUnderrun, 324
 - Status_Canceled, 324
 - Status_Canceling, 325
 - Status_ConnectedPipes, 325
 - Status_ControlNotOwned, 325

- Status_Crc, [325](#)
- Status_DataOverrun, [325](#)
- Status_DataToggleMismatch, [325](#)
- Status_DataUnderrun, [325](#)
- Status_DeviceLocked, [325](#)
- Status_DeviceNotFound, [325](#)
- Status_DeviceRemoved, [325](#)
- Status_DevNotResponding, [325](#)
- Status_EndpointHalted, [326](#)
- Status_ErrorBusy, [326](#)
- Status_ErrorShortTransfer, [326](#)
- Status_Fifo, [326](#)
- Status_FrameControlOwned, [326](#)
- Status_InternalHcError, [326](#)
- Status_InvalidDeviceHandle, [326](#)
- Status_InvalidHandle, [326](#)
- Status_InvalidParameter, [326](#)
- Status_InvalidPipeHandle, [326](#)
- Status_InvalidUrbFunction, [327](#)
- Status_IoPending, [327](#)
- Status_IoTimeout, [327](#)
- Status_IsochRequestFailed, [327](#)
- Status_LastUsbErrorMismatch, [327](#)
- Status_NoBandwidth, [327](#)
- Status_NoMemory, [327](#)
- Status_NoSuchDevice, [327](#)
- Status_NotAccessed, [327](#)
- Status_NotSupported, [327](#)
- Status_PidCheckFailure, [327](#)
- Status_PipeNotLinked, [328](#)
- Status_RequestFailed, [328](#)
- Status_RequestMutexFailed, [328](#)
- Status_RequestMutexTimeout, [328](#)
- Status_Stall, [328](#)
- Status_Unconfigured, [328](#)
- Status_UnexpectedPid, [328](#)
- ThrowCUsbExceptionNetOnError, [322](#)
- TxnGetSerialNumber, [322](#)
- TxnSetSerialNumber, [322](#)
- TxnTestMemoryReadAndCheck, [322](#)
- TxnTestMemoryWrite, [322](#)
- ValidKey, [322](#)
- WPAError_ScanningIsPending, [328](#)
- WriteEepromRegisterPreconfig, [323](#)
- WriteRegister, [323](#)
- WriteRegister32, [323](#)
- WriteRegisterArray, [323](#)
- WriteRegisterTimeSlot, [324](#)
- WriteRegisterValue, [324](#)
- CMcsUsbPointerContainer, [329](#)
- CMEA2100_256DacqGroupChannelSelectionNet, [329](#)
 - CMEA2100_256DacqGroupChannelSelectionNet, [329](#)
- CMEA2100x256FunctionNet, [329](#)
 - !CMEA2100x256FunctionNet, [330](#)
 - ~CMEA2100x256FunctionNet, [330](#)
 - CMEA2100x256FunctionNet, [330](#)
 - GetLayoutConfiguration, [330](#)
 - SetLayoutConfiguration, [330](#)
- CMeaAudioFunctionNet, [331](#)
 - CMeaAudioFunctionNet, [331](#)
 - GetAudioChannels, [332](#)
 - GetNumberOfAudioChannels, [332](#)
 - SetAudioChannels, [333](#)
- CMeaAudioFunctionNet::s_setaudionet, [693](#)
 - amplification, [693](#)
 - channel, [693](#)
- CMeaCleanDeviceNet, [334](#)
 - !CMeaCleanDeviceNet, [335](#)
 - ~CMeaCleanDeviceNet, [335](#)
 - CMeaCleanDeviceNet, [335](#)
 - GetCycle, [335](#)
 - GetCycles, [335](#)
 - GetMaxVoltage, [335](#)
 - GetMinVoltage, [336](#)
 - GetOutputVoltage, [336](#)
 - GetSlope, [336](#)
 - IsRunning, [336](#)
 - SetCycles, [336](#)
 - SetMaxVoltage, [337](#)
 - SetMinVoltage, [337](#)
 - SetSlope, [337](#)
 - Start, [337](#)
 - Stop, [337](#)
- CMeaCoatDeviceNet, [338](#)
 - !CMeaCoatDeviceNet, [339](#)
 - ~CMeaCoatDeviceNet, [339](#)
 - CMeaCoatDeviceNet, [339](#)
 - GetCurrentCycle, [339](#)
 - GetCycles, [339](#)
 - GetDuration, [340](#)
 - GetMaxCurrent, [340](#)
 - GetOffsetCurrent, [340](#)
 - GetOutputCurrent, [340](#)
 - GetPauseDuration, [340](#)
 - GetSlope, [341](#)
 - GetTimeInPause, [341](#)
 - GetTimeInPlateau, [341](#)
 - IsRunning, [341](#)
 - SetCycles, [341](#)
 - SetDuration, [342](#)
 - SetMaxCurrent, [342](#)
 - SetOffsetCurrent, [342](#)
 - SetPauseDuration, [342](#)
 - SetSlope, [343](#)
 - Start, [343](#)
 - Stop, [343](#)
- CMeaDeviceNet, [343](#)
 - ~CMeaDeviceNet, [345](#)
 - AnalogGain, [350](#)
 - CMeaDeviceNet, [345](#)
 - EnableChecksum, [345](#)
 - EnableDigitalIn, [346](#)
 - EnableTimestamp, [347](#)
 - Gain, [350](#)
 - GetAnalogGain, [347](#)

- GetEnumerationSpeed, 347
- GetGain, 347
- MeaAudioFunctionNet, 351
- MeaDigitalDataFunctionNet, 351
- MeaFeedbackFunctionNet, 351
- MeFunctionNet, 351
- SetDigitalOut, 347
- SetNumberOfAnalogChannels, 348
- SetNumberOfChannels, 349
- SetTriggerMaskValue, 349
- SetTriggerPeriod, 350
- W2100_FunctionNet, 351
- WClassicFunctionNet, 351
- CMeaDigitalDataFunctionNet, 351
 - CMeaDigitalDataFunctionNet, 352
 - GetDigitalData, 352
 - SetDigitalData, 352, 353
- CMeaFeedbackFunctionNet, 353
 - CMeaFeedbackFunctionNet, 354
 - FeedbackGetSampleTimerCount, 354
 - FeedbackSetAnalogSource, 354
 - FeedbackSetChannelFilter, 354
 - FeedbackSetDigitalMapping, 355
 - FeedbackSetFeedback, 355
 - FeedbackSetFilterOff, 355
 - FeedbackSetFilterParameter, 355
 - FeedbackSetFilterParameter32, 355
 - FeedbackSetGlobalChannelFilter, 355
 - FeedbackSetIIRFilterParameter, 355
 - FeedbackSetLogic, 355
 - FeedbackSetMkFilter, 356
 - FeedbackSetNumberOfLogics, 356
 - FeedbackSetNumberOfRateCounter, 356
 - FeedbackSetNumberOfRateDetectors, 356
 - FeedbackSetNumberOfSpikeDetectors, 356
 - FeedbackSetNumberOfTriggers, 356
 - FeedbackSetRateCounter, 356
 - FeedbackSetRateDetector, 356
 - FeedbackSetSpikeDetectorThreshold, 357
 - FeedbackSetTrigger, 357
- CMealImpedanceDeviceNet, 357
 - ~CMealImpedanceDeviceNet, 358
 - CMealImpedanceDeviceNet, 358
 - GetAdapterCode, 358
 - GetArraySize, 358
 - GetImpedanceTestFrequency, 358
 - GetReady, 358
 - GetResult, 358
 - SetImpedanceTestFrequency, 358
 - StartMeasurement, 358
- CMeasureTableDeviceNet, 359
 - CMeasureTableDeviceNet, 359
 - Sensor, 359
- CMeaSwitchDeviceNet, 360
 - ~CMeaSwitchDeviceNet, 360
 - CMeaSwitchDeviceNet, 360
 - GetNumber, 361
 - GetPattern, 361
 - GetPatternBool, 361
 - SetPattern, 361
 - SetPatternBool, 361
- CMeaUSBDeviceNet, 361
 - ~CMeaUSBDeviceNet, 362
 - CMeaUSBDeviceNet, 362
- CMeFunctionNet, 363
 - !CMeFunctionNet, 363
 - ~CMeFunctionNet, 363
 - CMeFunctionNet, 363
 - SetTransformer, 364
- CmosMea
 - CCMOSMeaDeviceNet, 110
- CmosMea
 - Mcs::Usb, 65
- CMOSMeaBathModeEnumNet
 - Mcs::Usb, 54
- CmosMeaHeadstage
 - Mcs::Usb, 64
- CMOSMeaHeadstage1NCBathCurrentEnumNet
 - Mcs::Usb, 54
- CMOSMeaHeadstage1NCCol2CurrentEnumNet
 - Mcs::Usb, 54
- CMOSMeaHeadstage1NChipTempEnumNet
 - Mcs::Usb, 55
- CMOSMeaHS1SidebandEnumNet
 - Mcs::Usb, 55
- CMOSMeaHS1TriggerStatusEnumNet
 - Mcs::Usb, 55
- CmosmealFB2
 - Mcs::Usb, 65
- CMOSMealFDigChannelEnumNet
 - Mcs::Usb, 55
- CMOSMealInterfaceADCEnumNet
 - Mcs::Usb, 56
- CmosMealInterfaceboard
 - Mcs::Usb, 64
- CMOSMeaPacketFrameContextGroupEnumNet
 - Mcs::Usb, 56
- CMOSMeaSTG1DACSignalEnumNet
 - Mcs::Usb, 56
- CMOSMeaValueUnitEnumNet
 - Mcs::Usb, 57
- CMultiBatteryChargerDeviceNet, 364
 - !CMultiBatteryChargerDeviceNet, 365
 - ~CMultiBatteryChargerDeviceNet, 365
 - CapacityTest, 365
 - ChannelReset, 366
 - CMultiBatteryChargerDeviceNet, 365
 - GetBatteryVoltage, 366
 - GetChannels, 366
 - GetChannelState, 366
 - GetChargeCapacity, 367
 - GetChargeCurrent, 367
 - GetChargingMode, 367
 - GetChargingPCoefficient, 368
 - GetDischargeCapacity, 368
 - GetDischargeCurrent, 368

- GetDischargeCurrentSetPoint, 368
- GetFinalDischargeVoltage, 369
- GetRatedCapacity, 369
- SetChargingMode, 369
- SetChargingPCoefficient, 370
- SetDischargeCurrentSetPoint, 370
- SetFinalDischargeVoltage, 370
- SetRatedCapacity, 370
- SetRatedCapacityVolatile, 372
- CMultiwellCallbackFunctionNet, 372
 - !CMultiwellCallbackFunctionNet, 373
 - ~CMultiwellCallbackFunctionNet, 373
 - CMultiwellCallbackFunctionNet, 373
 - GetPlateClampStateByHeadstage, 373
 - GetPlateClampStateByHeadstageEvent, 374
 - OnGetPlateClampStateByHeadstage, 374
- CMultiwellDeviceNet, 374
 - !CMultiwellDeviceNet, 376
 - ~CMultiwellDeviceNet, 376
 - ClosePlateClamp, 376
 - CMultiwellDeviceNet, 375
 - GetPlateClampLockState, 376
 - GetPlateClampState, 376
 - GetPlateMux, 377
 - GetPlateType, 377
 - GetPowerMuxPlate, 378
 - GetTouchPadEnable, 378
 - GetVolatileClampOffset, 378
 - IsPlateTypeValid, 379
 - LockPlateClamp, 379
 - OpenPlateClamp, 379
 - SetPlateMux, 379, 380
 - SetPlateType, 380
 - SetPowerMuxPlate, 380
 - SetTouchPadEnable, 381
 - SetVolatileClampOffset, 381
 - StopPlateClamp, 381
 - UnlockPlateClamp, 381
- CMultiwellOptoStimFunctionNet, 382
 - !CMultiwellOptoStimFunctionNet, 383
 - ~CMultiwellOptoStimFunctionNet, 383
 - CMultiwellOptoStimFunctionNet, 382
 - GetAbsMaxCurrentInMicroAmp, 383
 - GetColorRgb, 383
 - GetColorStr, 383
 - GetMaxDurationHighCurrentInMicroSec, 384
 - GetMaxDutyCycleHighCurrent, 384
 - GetPermanentCurrentInMicroAmp, 384
 - GetWaveLengthInNanometer, 384
 - SetAbsMaxCurrentInMicroAmp, 385
 - SetColorRgb, 385
 - SetColorStr, 385
 - SetMaxDurationHighCurrentInMicroSec, 385
 - SetMaxDutyCycleHighCurrent, 386
 - SetPermanentCurrentInMicroAmp, 386
 - SetWaveLengthInNanometer, 386
- CNF_GenDeviceNet, 386
 - ~CNF_GenDeviceNet, 387
 - CNF_GenDeviceNet, 387
 - Set_Values, 387
- COctoPotDeviceNet, 387
 - BurnAdcOffset, 388
 - BurnDacOffset, 389
 - COctoPotDeviceNet, 388
 - EnableChecksum, 389
 - EnableDigitalIn, 389
 - EnableTimestamp, 389
 - GetAdcOffset, 389
 - GetDacOffset, 389
 - PatternListStart, 389
 - RampStart, 389
 - ResetAdcOffset, 389
 - ResetDacOffset, 389
 - SetAdcOffset, 390
 - SetAmplificationSwitch, 390
 - SetBathclamp, 390
 - SetChannelSwitch, 390
 - SetDacAutoControl, 390
 - SetDacOffset, 390
 - SetDacValue, 390
 - SetNumberOfChannels, 390
 - SetOutputRate, 390
 - SetPatternListEntry, 391
 - SetPidParameter, 391
 - SetRampParameter, 391
 - SetSineParameter, 391
 - SineStart, 391
- CokuvisionStimulatorDeviceNet, 391
 - ~CokuvisionStimulatorDeviceNet, 392
 - CokuvisionStimulatorDeviceNet, 392
 - GetCheckVoltage, 392
 - GetCurrentFactor, 392
 - GetDACOffset, 393
 - GetMaxPower, 393
 - GetMaxVoltage, 393
 - GetPulseform, 393
 - GetRTC, 393
 - GetStimulatorStatus, 393
 - GetVoltage, 393
 - SetCheckVoltage, 394
 - SetCurrentFactor, 394
 - SetDACOffset, 394
 - SetMaxPower, 394
 - SetMaxVoltage, 394
 - SetPulseform, 394
 - SetRTC, 394
- Coldstart
 - CMcsUsbFactoryNet, 291
- CommaPositionA
 - CFilterCoefficientsNet::s_FilterAttributesNet, 693
 - Mcs::Usb, 63
- CommaPositionB
 - CFilterCoefficientsNet::s_FilterAttributesNet, 693
 - Mcs::Usb, 63
- CompareFirmware
 - CMcsUsbFactoryNet, 292

- CompareTo
 - HeadStageIDType, 678
- CompensateElectrodeOffset
 - CWarnerUssingFunctionNet, 615
- Connect
 - CMcsUsbNet, 312, 313
 - CRFFunctionNet, 439
- ConnectDevice
 - CRadioControlledDevicesNet, 434
- ConnectedImp
 - CPositionImpDeviceNet, 411
- ConnectImp
 - CPositionImpDeviceNet, 411
- ConnectSlave
 - CGilsonDeviceNet, 164
- ControlState
 - HeadStageIDTypeState, 681
- CornerFrequency
 - CFilterPropertyNet, 135
- CornerFrequencyHz
 - CFilterPropertyNet, 136
- Count
 - CMcsUsbListNet, 307
- CPatchServerDeviceNet, 395
 - CPatchServerDeviceNet, 395
 - Sensor, 395
- CPathIdentDeviceNet, 396
 - ~CPathIdentDeviceNet, 396
 - CPathIdentDeviceNet, 396
 - Measure, 396
 - Set_Values, 396
- CPedoterDeviceNet, 397
 - !CPedoterDeviceNet, 397
 - ~CPedoterDeviceNet, 397
 - CPedoterDeviceNet, 397
 - GetCommand, 397
 - SetCommand, 398
- CPeristalticPumpDeviceNet, 398
 - ~CPeristalticPumpDeviceNet, 399
 - CPeristalticPumpDeviceNet, 399
 - McsBus_MotorControl, 399
- CPgaDeviceNet, 399
 - ~CPgaDeviceNet, 400
 - ApplyGains, 400
 - CPgaDeviceNet, 400
 - DefineAmplification, 400
 - DefineFrequencyRange, 400
 - DefineNumAmplifications, 400
 - DefineNumFrequencyRanges, 400
 - GetAmplification, 400
 - GetFrequencyRange, 401
 - GetGain, 401
 - GetNumAmplifications, 401
 - GetNumFrequencyRanges, 401
 - SetGain, 401
- CPositionIIDeviceNet, 401
 - !CPositionIIDeviceNet, 403
 - ~CPositionIIDeviceNet, 403
 - CPositionIIDeviceNet, 403
 - GetCoilCommunication, 403
 - GetDebugData, 403
 - GetEventData, 404
 - GetImplantCurrentSetpoint, 404
 - GetImplantResult, 405
 - GetImplantState, 405
 - GetOnOff, 405
 - GetPowerStrength, 406
 - GetRTC, 406
 - GetStateDebugData, 406
 - GetStateEventData, 407
 - RFFunction, 410
 - SetImplantCurrentSetpoint, 407
 - SetPowerStrength, 407
 - SetRTC, 407
 - SetStateDebugData, 409
 - SetStateEventData, 409
 - SwitchOnOff, 409
- CPositionImpDeviceNet, 410
 - !CPositionImpDeviceNet, 411
 - ~CPositionImpDeviceNet, 411
 - ConnectedImp, 411
 - ConnectImp, 411
 - CPositionImpDeviceNet, 411
 - GetDeviceList, 411
 - GetImpld, 412
 - GetRFFrequency, 412
 - SetDeviceList, 412
 - SetImpld, 412
 - SetRFFrequency, 413
- CPPCDeviceNet, 413
 - CPPCDeviceNet, 413
 - McsBus, 414
 - McsBus_MotorControl, 414
 - McsBus_Sensor, 414
 - PPCFunction, 414
- CPPPCFunctionNet, 414
 - !CPPPCFunctionNet, 416
 - ~CPPPCFunctionNet, 416
 - CPPPCFunctionNet, 415
 - FirePressurePulse, 416
 - GetAnalogVoltage, 416
 - GetAnalogVoltageRange, 416
 - GetDigitalIn, 418
 - GetPressureRange, 418
 - GetPumpModeType, 418
 - GetPumpSpeedUnit, 419
 - GetSupplyVoltage, 419
 - GetValveActive, 419
 - IsBusy, 419
 - LoadPressure, 421
 - MeasureReservoir, 421
 - SetAnalogVoltageRange, 421
 - SetPressureOffset, 421
 - SetPressureRange, 421
 - SetPumpModeType, 422
 - SetPumpSpeedUnit, 422

- SetValveActive, [422](#)
- CPPS_DeviceNet, [423](#)
 - CPPS_DeviceNet, [423](#)
 - McsBus, [423](#)
 - McsBus_MotorControl, [423](#)
 - McsBus_Sensor, [423](#)
 - PPS_Function, [423](#)
- CPPS_FunctionNet, [424](#)
 - CPPS_FunctionNet, [424](#), [425](#)
 - GetAnalogVoltage, [425](#)
 - GetAnalogVoltages, [425](#)
 - GetBubbleState, [425](#)
 - GetDigitalIn, [425](#)
 - GetPumpCouple, [425](#)
 - GetPumpEnableSpeedRatio, [425](#)
 - GetPumpFastOnOff, [425](#)
 - GetPumpFastSpeed, [425](#)
 - GetPumpFunctionSpeeds, [426](#)
 - GetPumpManualOnOff, [426](#)
 - GetPumpMaxSpeed, [426](#)
 - GetPumpModeType, [426](#)
 - GetPumpSpeedRatio, [426](#)
 - GetPumpSpeedUnit, [426](#)
 - GetSupplyVoltage, [426](#)
 - GetUseBubble, [426](#)
 - SetAnalogVoltages, [426](#)
 - SetPumpCouple, [427](#)
 - SetPumpEnableSpeedRatio, [427](#)
 - SetPumpFastOnOff, [427](#)
 - SetPumpFastSpeed, [427](#)
 - SetPumpFunctionSpeeds, [427](#)
 - SetPumpManualOnOff, [427](#)
 - SetPumpMaxSpeed, [427](#)
 - SetPumpModeType, [427](#)
 - SetPumpSpeedRatio, [428](#)
 - SetPumpSpeedUnit, [428](#)
 - SetUseBubble, [428](#)
- CPPSDeviceNet, [428](#)
 - CPPSDeviceNet, [429](#)
- CProgramPressureCurveNet, [429](#)
 - !CProgramPressureCurveNet, [430](#)
 - ~CProgramPressureCurveNet, [429](#)
 - CProgramPressureCurveNet, [429](#)
 - GetRepeats, [430](#)
 - Program, [430](#)
 - SetRepeats, [430](#)
- CPulseGeneratorFunctionNet, [430](#)
 - !CPulseGeneratorFunctionNet, [431](#)
 - ~CPulseGeneratorFunctionNet, [431](#)
 - CPulseGeneratorFunctionNet, [431](#)
 - GetModeSelect, [432](#)
 - GetPeriod, [432](#)
 - GetPulseLength, [432](#)
 - SetModeSelect, [433](#)
 - SetPeriod, [433](#)
 - SetPulseLength, [433](#)
- CRadioControlledDevicesNet, [433](#)
 - ConnectDevice, [434](#)
 - CRadioControlledDevicesNet, [434](#)
 - DisConnectDevice, [434](#)
 - GetDeviceNames, [434](#)
 - GetFrequency, [435](#)
 - HasRadioControl, [435](#)
 - SetFrequency, [435](#)
 - StillConnected, [435](#)
- CreateSideband
 - CStimulusFunctionNet, [559](#)
- CreateWirelessHeadstageSerialNumberString
 - CWirelessBaseFunctionNet, [663](#)
- CRegionOfInterestRect
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [435](#)
- CRetinaLedDeviceNet, [436](#)
 - ~CRetinaLedDeviceNet, [437](#)
 - AddLoopEntry, [437](#)
 - AddTableEntry, [437](#)
 - ClearTable, [437](#)
 - CRetinaLedDeviceNet, [437](#)
 - GetTablepointer, [437](#)
 - SetLED, [437](#)
 - SetLumi, [437](#)
 - SetPersistency, [437](#)
 - SetRepeat, [438](#)
 - SetTablepointer, [438](#)
 - SetTrigger, [438](#)
- CRFFFunctionNet, [438](#)
 - !CRFFFunctionNet, [439](#)
 - ~CRFFFunctionNet, [439](#)
 - Connect, [439](#)
 - CRFFFunctionNet, [439](#)
 - GetAvailableDeviceList, [440](#)
 - GetAvailableDeviceListEx, [440](#)
 - GetAvailableStateList, [440](#)
 - GetAvailableStateListEx, [440](#)
 - GetBaseFrequency, [441](#)
 - GetConnectedDevice, [441](#)
 - GetState, [441](#)
 - GetTestMode, [441](#)
 - GetWorkingFrequency, [442](#)
 - SetBaseFrequency, [442](#)
 - SetTestMode, [442](#)
 - SetWorkingFrequency, [442](#)
- CRobo_FYIProgram_FunctionNet, [443](#)
 - CRobo_FYIProgram_FunctionNet, [443](#)
 - GetLength, [443](#)
 - GetState, [444](#)
 - GetValve1, [444](#)
 - GetValve2, [444](#)
 - SetLength, [444](#)
 - SetValve1, [444](#)
 - SetValve2, [444](#)
 - Start, [444](#)
- CRobo_FYITemp_FunctionNet, [444](#)
 - CRobo_FYITemp_FunctionNet, [445](#)
 - GetICoeff, [445](#)
 - GetMaxPower, [445](#)

- GetPCoeff, 445
- GetRegulatorOnOff, 445
- GetSollTemp, 445
- SetICoeff, 446
- SetMaxPower, 446
- SetPCoeff, 446
- SetRegulatorOnOff, 446
- SetSollTemp, 446
- CRoboDacqNet, 446
 - CancelTableLoop, 449
 - CancelTableLoopAndStopTable, 449
 - ClampAmpRestart, 449
 - CRoboDacqNet, 449
 - DoRamp, 449
 - Emu_GetCellCapacity, 450
 - Emu_GetCellPotential, 450
 - Emu_GetCellResists, 450
 - Emu_GetElectrodeResists, 450
 - Emu_GetNoise, 450
 - Emu_SetCellCapacity, 450
 - Emu_SetCellPotential, 450
 - Emu_SetCellResists, 450
 - Emu_SetElectrodeResists, 450
 - Emu_SetNoise, 450
 - GetAllDigout, 451
 - GetCapacityC, 451
 - GetCapacityV, 451
 - GetCapacityX, 451
 - GetClampAmpSerialNumber, 451
 - GetCommand, 451
 - GetConfigurationBit, 451
 - GetConfigurationBitAxc, 451
 - GetConfigurationBitBlu_Led, 451
 - GetConfigurationBitBlu_LedToggleFast, 451
 - GetConfigurationBitBlu_LedToggleSlow, 451
 - GetConfigurationBitCC_Gen, 452
 - GetConfigurationBitCV_Gen, 452
 - GetConfigurationBitRC_Gen, 452
 - GetConfigurationBitRed_Led, 452
 - GetConfigurationBitRed_LedSaturation, 452
 - GetConfigurationBitRed_LedToggleFast, 452
 - GetConfigurationBitRed_LedToggleSlow, 452
 - GetConfigurationBitRelais, 452
 - GetConfigurationBitRV_Gen, 452
 - GetConfigurationBits, 452
 - GetConfigurationBitStream, 452
 - GetConfigurationBitSupply, 453
 - GetCrossTalkOffset, 453
 - GetCrossTalkOptimum, 453
 - GetDigout, 453
 - GetDisplayText, 453
 - GetDownsampleFactor, 453
 - GetFilter, 453
 - GetFilterCoeffs, 453
 - GetIC, 453
 - GetIClamp, 453
 - GetICOffset, 453
 - GetIGain, 454
 - GetNIC_MS, 454
 - GetNUC_MS, 454
 - GetNUV_MS, 454
 - GetPGain, 454
 - GetRecordingNumber, 454
 - GetResistanceC, 454
 - GetResistanceV, 454
 - GetScreen, 454
 - GetSimulation, 454
 - GetUC, 454
 - GetUClamp, 455
 - GetUCOffset, 455
 - GetUpdateDisplay, 455
 - GetUV, 455
 - GetUVOffset, 455
 - GetXGain, 455
 - RunTable, 455
 - SetAllDigout, 455
 - SetCommand, 455
 - SetConfigurationBit, 455
 - SetConfigurationBitAxc, 456
 - SetConfigurationBitBlu_Led, 456
 - SetConfigurationBitBlu_LedToggleFast, 456
 - SetConfigurationBitBlu_LedToggleSlow, 456
 - SetConfigurationBitCC_Gen, 456
 - SetConfigurationBitCV_Gen, 456
 - SetConfigurationBitRC_Gen, 456
 - SetConfigurationBitRed_Led, 456
 - SetConfigurationBitRed_LedSaturation, 456
 - SetConfigurationBitRed_LedToggleFast, 457
 - SetConfigurationBitRed_LedToggleSlow, 457
 - SetConfigurationBitRelais, 457
 - SetConfigurationBitRV_Gen, 457
 - SetConfigurationBitStream, 457
 - SetConfigurationBitSupply, 457
 - SetCrossTalkOffset, 457
 - SetCrossTalkOptimum, 457
 - SetDigout, 457
 - SetDisplayText, 458
 - SetDownsampleFactor, 458
 - SetFilter, 458
 - SetFilterCoeffs, 458
 - SetIClamp, 458
 - SetICOffset, 458
 - SetIGain, 458
 - SetNoFilterCoeffs, 458
 - SetPGain, 458
 - SetRecordingNumber, 459
 - SetScreen, 459
 - SetSimulation, 459
 - SetTriggerMaskValue, 459
 - SetUClamp, 459
 - SetUCOffset, 459
 - SetUVOffset, 459
 - SetXGain, 459
 - StopTable, 459, 460
 - Table_Wait, 460
 - TableDefBegin, 460

- TableDefEnd, [460](#)
- TriggerMask_Default, [460](#)
- TriggerValue_MoveAbs, [460](#)
- TriggerValue_StartQueue, [460](#)
- UpdateDisplay, [460](#)
- VirtualDevice_ContinuousDacq, [460](#)
- VirtualDevice_TableRun, [460](#)
- CRoboDeviceNet, [461](#)
 - ~CRoboDeviceNet, [464](#)
 - Axes_I, [472](#)
 - Axes_X, [472](#)
 - Axes_Y, [472](#)
 - Axes_Z, [473](#)
 - Axis_I, [473](#)
 - Axis_X, [473](#)
 - Axis_Y, [473](#)
 - Axis_Z, [473](#)
 - CancelPoolLoop, [464](#)
 - CancelPoolLoopAndStopMovement, [464](#)
 - CRoboDeviceNet, [464](#)
 - EnableQueue, [465](#)
 - FindReference, [465](#)
 - GetAirpressure, [465](#)
 - GetAirpressureLimit, [465](#)
 - GetAirValve, [465](#)
 - GetCurrentAirvalve, [465](#)
 - GetCurrentAirvalveLimit, [465](#)
 - GetCurrentPosition, [466](#)
 - GetErrorAirpressure, [466](#)
 - GetErrorCurrentAirvalve, [466](#)
 - GetErrorVoltage12V, [466](#)
 - GetErrorVoltage5V, [466](#)
 - GetErrorVoltageAirvalve, [466](#)
 - GetErrorVoltageRs485A, [466](#)
 - GetErrorVoltageRs485B, [466](#)
 - GetErrorVoltageValves, [467](#)
 - GetInMovement, [467](#)
 - GetMinPressure, [467](#)
 - GetMovementError, [467](#)
 - GetVoltage12V, [467](#)
 - GetVoltage12VLimit, [467](#)
 - GetVoltage5V, [467](#)
 - GetVoltage5VLimit, [467](#)
 - GetVoltageAirvalve, [467](#)
 - GetVoltageAirvalveLimit, [467](#)
 - GetVoltageRs485A, [468](#)
 - GetVoltageRs485ALimit, [468](#)
 - GetVoltageRs485B, [468](#)
 - GetVoltageRs485BLimit, [468](#)
 - GetVoltageValves, [468](#)
 - GetVoltageValvesLimit, [468](#)
 - IsQueueEnabled, [468](#)
 - IsQueueStarted, [468](#)
 - McsBus, [476](#)
 - McsBus_MotorControl, [476](#)
 - McsBus_XY, [473](#)
 - McsBus_ZI, [473](#)
 - MoveAbs, [468](#), [469](#)
 - NullCommand, [470](#)
 - RoboError_AnotherMaster, [473](#)
 - RoboError_Base, [474](#)
 - RoboError_CannotEscapeEndSwitch, [474](#)
 - RoboError_CommandAlreadyInProgress, [474](#)
 - RoboError_CommandNotPossible, [474](#)
 - RoboError_CommunicationTimeout, [474](#)
 - RoboError_DacqNotReady, [474](#)
 - RoboError_DLLMovementTimeout, [474](#)
 - RoboError_FindReferenceMethod, [474](#)
 - RoboError_GilsonCommandPending, [474](#)
 - RoboError_GilsonTimeout, [474](#)
 - RoboError_GilsonWrondID, [475](#)
 - RoboError_McsBus_UnknownCommand, [475](#)
 - RoboError_NoEndSwitch, [475](#)
 - RoboError_NoMoreData, [475](#)
 - RoboError_NoReference, [475](#)
 - RoboError_NoSpeedOrAcceleration, [475](#)
 - RoboError_OverPressure, [475](#)
 - RoboError_ParameterNotAllowed, [475](#)
 - RoboError_PeristalticTimeout, [475](#)
 - RoboError_Phase0OutOfRange, [475](#)
 - RoboError_PollLoopCanceled, [476](#)
 - RoboError_PollLoopCanceledAndStopMovement, [476](#)
 - RoboError_Pressure, [476](#)
 - RoboError_RangeExceeded, [476](#)
 - RoboError_StateChangeNotPossible, [476](#)
 - RoboError_Timeout, [476](#)
 - RoboError_UnknownCommand, [476](#)
 - RoboMainLowLevelCommand, [477](#)
 - RoboStatusEvent, [477](#)
 - SetAirpressureLimit, [470](#)
 - SetAirValve, [470](#)
 - SetCurrentAirvalveLimit, [470](#)
 - SetCurrentAndAir, [470](#)
 - SetInMovement, [470](#)
 - SetMinPressure, [471](#)
 - SetVoltage12VLimit, [471](#)
 - SetVoltage5VLimit, [471](#)
 - SetVoltageAirvalveLimit, [471](#)
 - SetVoltageRs485ALimit, [471](#)
 - SetVoltageRs485BLimit, [471](#)
 - SetVoltageValvesLimit, [471](#)
 - StartQueue, [471](#)
 - StopMovement, [472](#)
 - WaitTimer, [472](#)
- CRoboDeviceNet::RoboMainLowLevelCommands, [685](#)
 - FindReferencePhase0, [686](#)
 - GetAxisConfig, [686](#)
 - GetHWConfig, [686](#)
 - GetHWRevision, [686](#)
 - GetMaxNoPressure, [686](#)
 - GetMaxNoPressureWaitTime, [686](#)
 - GetMaxPressureWaitTime, [686](#)
 - GetMinNoPressureWaitTime, [686](#)
 - GetMinPressure, [686](#)
 - GetMinPressureWaitTime, [687](#)

- GetParameter, [687](#)
- GetPhases, [687](#)
- GetSearchReferenceFastAccel, [687](#)
- GetSearchReferenceFastSpeed, [687](#)
- GetSearchReferenceFineAccel, [687](#)
- GetSearchReferenceFineSpeed, [687](#)
- GetSearchReferenceMethod, [688](#)
- GetSearchReferenceMoveOut, [688](#)
- GetSearchReferenceOffsetPos, [688](#)
- GetUserParameter, [688](#)
- HasRef, [689](#)
- SetAxisConfig, [689](#)
- SetHWConfig, [689](#)
- SetHWRRevision, [689](#)
- SetMaxNoPressure, [689](#)
- SetMaxNoPressureWaitTime, [689](#)
- SetMaxPressureWaitTime, [689](#)
- SetMinNoPressureWaitTime, [689](#)
- SetMinPressure, [689](#)
- SetMinPressureWaitTime, [689](#)
- SetParameter, [690](#)
- SetSearchReferenceFastAccel, [690](#)
- SetSearchReferenceFastSpeed, [690](#)
- SetSearchReferenceFineAccel, [690](#)
- SetSearchReferenceFineSpeed, [690](#)
- SetSearchReferenceMethod, [690](#)
- SetSearchReferenceMoveOut, [691](#)
- SetSearchReferenceOffsetPos, [691](#)
- SetUserParameter, [691](#)
- CRoboFluidDeviceNet, [477](#)
 - ~CRoboFluidDeviceNet, [478](#)
 - CloseAllValves, [478](#)
 - CRoboFluidDeviceNet, [478](#)
 - GetPumpSpeed, [478](#)
 - GetSingleValve, [478](#)
 - GetValve, [478](#)
 - IsPumpMotorOn, [479](#)
 - m_pMcsBus_MotorControlNet, [480](#)
 - m_pRoboFluidDevice, [480](#)
 - McsBus_MotorControl, [480](#)
 - PumpOff, [479](#)
 - PumpOn, [479](#)
 - SetPumpSpeed, [479](#)
 - SetSingleValve, [479](#)
 - SetValve, [479](#)
- CRoboInjectDeviceNet, [480](#)
 - CRoboInjectDeviceNet, [481](#)
- CRoboocyte2DeviceNet, [481](#)
 - CRoboocyte2DeviceNet, [482](#)
 - GetAxisLED, [482](#)
 - GetGilsonDevice, [482](#)
 - GetMcsBus_Extension, [482](#)
 - GetRoboDacq, [482](#)
 - GetRoboFluidDevice, [482](#)
 - SetAxisLED, [482](#)
- CRoboStatorDeviceNet, [482](#)
 - CRoboStatorDeviceNet, [484](#)
 - FindReferenceI, [484](#)
 - FindReferenceXY, [484](#)
 - FindReferenceZ, [484](#)
 - GetCurrentPositionI, [484](#)
 - GetCurrentPositionXY, [484](#)
 - GetCurrentPositionZ, [485](#)
 - HasRefI, [485](#)
 - HasRefXY, [485](#)
 - HasRefZ, [485](#)
 - MoveAbsI, [485](#)
 - MoveAbsXY, [485](#)
 - MoveAbsZ, [485, 486](#)
 - RoboMainStatorLowLevelCommand, [488](#)
 - SetAccelerationI, [486](#)
 - SetAccelerationNativeI, [486](#)
 - SetAccelerationNativeXY, [486](#)
 - SetAccelerationNativeZ, [486](#)
 - SetAccelerationXY, [486](#)
 - SetAccelerationZ, [486](#)
 - SetCurrentAndAirXY, [486](#)
 - SetSpeedI, [487](#)
 - SetSpeedNativeI, [487](#)
 - SetSpeedNativeXY, [487](#)
 - SetSpeedNativeZ, [487](#)
 - SetSpeedXY, [487](#)
 - SetSpeedZ, [487](#)
 - SetVelocityI, [487](#)
 - SetVelocityXY, [487](#)
 - SetVelocityZ, [487](#)
 - StopMovementI, [488](#)
 - StopMovementXY, [488](#)
 - StopMovementZ, [488](#)
- CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands, [692](#)
 - FindReferencePhase0XY, [692](#)
- CSafeISDeviceNet, [488](#)
 - ~CSafeISDeviceNet, [489](#)
 - CSafeISDeviceNet, [489](#)
 - DacqDevice, [491](#)
 - FluidControlDevice, [491](#)
 - RoboDevice, [491](#)
 - SetAdcChannels, [490](#)
 - SetAdcSamplePos, [490](#)
 - SetDacMode, [490](#)
 - SetDacPeriode, [490](#)
 - SetDacPulseform, [490](#)
 - SetSwitches, [491](#)
- csCapacityTestDischarge
 - Mcs::Usb, [67](#)
- csCapacityTestPrecharge
 - Mcs::Usb, [67](#)
- csCharge
 - Mcs::Usb, [67](#)
- CSCUDacqGroupChannelSelectionNet, [491](#)
 - CSCUDacqGroupChannelSelectionNet, [492](#)
- CSCUFunctionNet, [492](#)
 - !CSCUFunctionNet, [495](#)
 - ~CSCUFunctionNet, [495](#)
 - AutomaticAnalogOut, [495](#)

- CSCUFunctionNet, [494](#)
- EnableAnalogOut, [495](#)
- GetAnalogOutADCRange, [495](#)
- GetAnalogOutChannels, [495](#)
- GetAnalogOutDACRange, [496](#)
- GetAvailableHeadstages, [496](#)
- GetAvailableHeadstagesEvent, [506](#)
- GetFilterProperties, [496](#)
- GetFilterProperty, [497](#)
- GetHeadstageAdcBits, [497](#)
- GetHeadstageAdcRangeInMicroVolt, [497](#)
- GetHeadstageDacBits, [498](#)
- GetHeadstageDacCurrentRangeInMicroAmpere, [498](#)
- GetHeadstageDacCurrentResolutionInNanoAmpere, [498](#)
- GetHeadstageDacVoltageRangeInMilliVolt, [499](#)
- GetHeadstageDacVoltageResolutionInMicroVolt, [499](#)
- GetHeadstageGainInPer mille, [499](#)
- GetHeadstageID, [500](#)
- GetHeadstageNumberOfAnalogChannels, [500](#)
- GetHeadstageNumberOfStimulationChannels, [500](#)
- GetHeadstagePowerStateAtStart, [501](#)
- GetHeadstageSamplerate, [501](#)
- GetHeadstageSerialNumber, [501](#)
- GetMaxNumberOfHeadstages, [502](#)
- GetMaxStimulusChannelsPerHeadstage, [502](#)
- GetReferenceElectrodeMode, [502](#)
- GetReferenceElectrodeSwitchState, [502](#)
- HasAnalogOut, [503](#)
- HasGalvanicIsolation, [503](#)
- HasHSPowerSwitch, [503](#)
- IsAnalogOutEnabled, [503](#)
- IsAutomaticAnalogOut, [503](#)
- IsHeadstageAvailable, [503](#)
- IsHeadstageAvailableEvent, [507](#)
- IsHSPowered, [504](#)
- IsInDacqLegacyMode, [504](#)
- OnGetAvailableHeadstages, [504](#)
- OnIsHeadstageAvailable, [504](#)
- PowerHS, [504](#)
- SetAnalogOutADCRange, [505](#)
- SetAnalogOutChannels, [505](#)
- SetAnalogOutDACRange, [505](#)
- SetDacqLegacyMode, [505](#)
- SetHeadstagePowerStateAtStart, [506](#)
- SetReferenceElectrodeMode, [506](#)
- SetReferenceElectrodeSwitchState, [506](#)
- csDischarge
 - Mcs::Usb, [67](#)
- CSerialPortNet, [507](#)
 - CSerialPortNet, [507](#)
 - GetBytesAvailable, [507](#)
 - Receive, [508](#)
 - ReceiveString, [508](#)
 - Send, [508](#)
- csError
 - Mcs::Usb, [67](#)
- csIdleChargeFinished
 - Mcs::Usb, [67](#)
- csIdleNoBattery
 - Mcs::Usb, [67](#)
- csRefreshBattery
 - Mcs::Usb, [67](#)
- CStg200xBasicNet, [508](#)
 - ~CStg200xBasicNet, [513](#)
 - GetAnalogRanges, [513](#)
 - GetAnalogResolution, [513](#)
 - GetAutocalibrationDisabled, [513](#)
 - GetAvailableMemory, [515](#)
 - GetBlankingEnable, [515](#)
 - GetCurrentRangeByIndex, [516](#)
 - GetCurrentRangeInNanoAmp, [516](#)
 - GetCurrentResolutionInNanoAmp, [516](#)
 - GetDacAmplificationFactor, [517](#)
 - GetDACResolution, [517](#)
 - GetDiginValue, [517](#)
 - GetDigoutMode, [517](#)
 - GetDigoutValue, [517](#)
 - GetElectrodeDacMux, [518](#)
 - GetElectrodeEnable, [518](#), [519](#)
 - GetElectrodeMode, [519](#), [520](#)
 - GetEnableAmplifierProtectionSwitch, [520](#)
 - GetExternalElectrodeEnable, [521](#)
 - GetFAAmplification, [521](#)
 - GetHeadstage, [521](#)
 - GetListModelIndexRange, [521](#)
 - GetListModelTriggerSource, [522](#)
 - GetNumberOfAnalogChannels, [522](#)
 - GetNumberOfHWDACPaths, [522](#)
 - GetNumberOfStimulationElectrodes, [523](#)
 - GetNumberOfStimulationSourcesPerElectrode, [523](#)
 - GetNumberOfSyncoutChannels, [523](#)
 - GetNumberOfTriggerInputs, [523](#)
 - GetOutputRate, [523](#)
 - GetStgProgramInfo, [523](#), [524](#)
 - GetStgVersionInfo, [524](#)
 - GetSyncoutMap, [524](#)
 - GetTotalMemory, [525](#)
 - GetTriggerSource, [525](#)
 - GetVoltageRangeInMicroVolt, [525](#)
 - GetVoltageResolutionInMicroVolt, [525](#)
 - ListModeSendStart, [526](#)
 - ListModeSendStop, [526](#)
 - SendStart, [526](#)
 - SendStop, [526](#)
 - SetAutocalibrationDisabled, [527](#)
 - SetBlankingEnable, [527](#), [528](#)
 - SetCurrentMode, [528](#)
 - SetCurrentRangeByIndex, [529](#)
 - SetDacAmplificationFactor, [529](#)
 - SetDigoutMode, [529](#)
 - SetDigoutValue, [529](#)
 - SetElectrodeDacMux, [530](#), [531](#)

- SetElectrodeEnable, [532](#), [533](#)
- SetElectrodeMode, [534](#), [535](#)
- SetEnableAmplifierProtectionSwitch, [536](#), [537](#)
- SetExternalElectrodeEnable, [537](#), [538](#)
- SetFAAmplification, [538](#)
- SetHeadstage, [538](#)
- SetListModelIndexRange, [538](#)
- SetListModelTriggerSource, [539](#)
- SetMeasurementMode, [540](#)
- SetOutputRate, [540](#)
- SetStgProgramInfo, [540](#)
- SetSyncoutMap, [541](#)
- SetTriggerSource, [541](#)
- SetVoltageMode, [541](#)
- CStg200xDownloadBasicNet, [542](#)
 - ClearChannelData, [543](#)
 - ClearSyncData, [544](#)
 - DisableAutoReset, [544](#)
 - EnableAutoReset, [544](#)
 - ForceStatusEvent, [544](#)
 - GetMemoryUsageDAC, [544](#)
 - GetMemoryUsageSyncout, [544](#)
 - GetSweepCount, [545](#)
 - GetTrigger, [545](#)
 - ResetStatus, [546](#)
 - SendChannelData, [546](#)
 - SendSyncData, [546](#)
 - SetupRetriggerMode, [547](#)
 - SetupTrigger, [548](#)
 - SetupTriggerSingle, [548](#)
 - Stimulus, [549](#)
- CStg200xDownloadNet, [549](#)
 - ~CStg200xDownloadNet, [551](#)
 - ClearChannel_PrepareAndSendData, [551](#)
 - CStg200xDownloadNet, [550](#)
 - DisableMultiFileMode, [552](#)
 - EnableMultiFileMode, [552](#)
 - GetModuleCurrent, [552](#)
 - GetModuleTemp, [552](#)
 - MwPollStatusEvent, [555](#)
 - PrepareAndAppendData, [552](#)
 - PrepareAndSendData, [553](#)
 - QueryTriggerstatus, [554](#)
 - SendSegmentDefine, [554](#)
 - SendSegmentSelect, [554](#)
 - SendSegmentStart, [555](#)
 - SetOutputMap, [555](#)
 - Stg200xPollStatusEvent, [556](#)
- CStimulusFunctionNet, [556](#)
 - ClearChannel_PrepareAndSendData, [558](#)
 - ClearChannelData, [558](#)
 - ClearMultiplexedData, [558](#)
 - ClearSyncData, [558](#)
 - CreateSideband, [559](#)
 - CStimulusFunctionNet, [557](#), [558](#)
 - ForceStatusEvent, [559](#)
 - GetAvailableMemory, [560](#)
 - GetCurrentRangeInNanoAmp, [560](#)
 - GetCurrentResolutionInNanoAmp, [560](#)
 - GetDACResolution, [560](#)
 - GetMultiplexedDataChannelsInBlock, [561](#)
 - GetNumberOfAnalogChannels, [561](#)
 - GetTotalMemory, [561](#)
 - GetVoltageRangeInMicroVolt, [561](#)
 - GetVoltageResolutionInMicroVolt, [562](#)
 - PollStatusEvent, [567](#)
 - PrepareAndAppendData, [562](#)
 - PrepareAndSendData, [563](#)
 - PrepareData, [564](#)
 - SendMultiplexedData, [564](#)
 - SendPreparedData, [564](#)
 - SendStart, [564](#)
 - SendStop, [565](#)
 - SetupTrigger, [565](#)
 - SetupTriggerSingle, [566](#)
 - StartPoll, [566](#)
 - StopPoll, [566](#)
- CStimulusFunctionNet::SidebandData, [694](#)
 - !SidebandData, [694](#)
 - ~SidebandData, [694](#)
 - Duration, [694](#)
 - Sideband, [694](#)
 - SidebandData, [694](#)
- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [695](#)
 - !StimulusDeviceDataAndUnrolledData, [696](#)
 - ~StimulusDeviceDataAndUnrolledData, [696](#)
 - DeviceData, [696](#)
 - DeviceDataLength, [696](#)
 - StimulusDeviceDataAndUnrolledData, [696](#)
 - UnrolledAmplitude, [696](#)
 - UnrolledDuration, [696](#)
 - UnrolledSync, [696](#)
- Csw2to64DeviceNet, [567](#)
 - ~Csw2to64DeviceNet, [568](#)
 - Csw2to64DeviceNet, [568](#)
 - GetChannel, [568](#)
 - GetChannels, [568](#)
 - GetNumber, [568](#)
 - SetChannel, [568](#)
 - SetChannels, [569](#)
- CTcxDeviceNet, [569](#)
 - ~CTcxDeviceNet, [572](#)
 - CalibrateThermocouple, [572](#)
 - CTcxDeviceNet, [571](#)
 - FactoryReset, [572](#)
 - GetBoardTemp, [572](#)
 - GetCalibration, [572](#)
 - GetCalibrationDecp, [572](#)
 - GetCalibrationMax, [572](#)
 - GetCalibrationMin, [572](#)
 - GetCurrent, [572](#)
 - GetD, [573](#)
 - GetDDecp, [573](#)
 - GetDevice, [573](#)
 - GetDeviceType, [573](#)

GetDevname, 573
GetDMax, 573
GetDMin, 573
GetDuty, 573
GetEnableHeaterLimit, 574
GetEnableThermocouple, 574
GetHasThermocouple, 574
GetHeaterLimit, 574
GetHeaterTemp, 574
GetI, 574
GetIDecp, 574
GetIMax, 575
GetIMin, 575
GetIOut, 575
GetMaxHeaterPowerMultiwell, 575
GetMaxP, 575
GetMaxpDecp, 575
GetMaxpMax, 575
GetMaxpMin, 576
GetNumControlChannels, 576
GetNumDevices, 576
GetNumMeasureChannels, 576
GetOnOff, 576
GetP, 576
GetPDecp, 576
GetPMax, 576
GetPMin, 576
GetPOut, 576
GetPwrOut, 577
GetPwrSet, 577
GetRes1, 577
GetRes2, 577
GetResS, 577
GetResX, 577
GetROut, 577
GetSensorType, 578
GetSetpoint, 578
GetSetpointDecp, 578
GetSetpointMax, 578
GetSetpointMin, 578
GetThermocoupleCalibration, 578
GetThermocoupleNanovoltPerKelvin, 578
GetThermocoupleReferenceTemp, 579
GetThermocoupleTemp, 579
GetThermocoupleTempAbs, 579
GetUnit, 579
GetUOut, 579
GetValue, 579
GetValueHires, 580
GetVolti, 580
SetCalibration, 580
SetD, 580
SetDevice, 580
SetDeviceType, 580
SetDevname, 580
SetEnableHeaterLimit, 580
SetEnableThermocouple, 581
SetHeaterLimit, 581
SetI, 581
SetMaxHeaterPowerMultiwell, 581
SetMaxP, 581
SetOnOff, 581
SetP, 582
SetSensorType, 582
SetSetpoint, 582
SetThermocoupleNanovoltPerKelvin, 582
CTEERFunctionNet, 582
 !CTEERFunctionNet, 585
 ~CTEERFunctionNet, 585
 CancelInternalCalibration, 585
 CTEERFunctionNet, 584
 GetAdapterCode, 585
 GetAdcOffsetU1, 585
 GetAdcOffsetU2, 585
 GetAmplitude_nA, 585
 GetBytesPerSample, 586
 GetClampMode, 586
 GetControllerParams, 586
 GetCurrentEnable, 586
 GetDacZero, 587
 GetFrameErrorCounter, 587
 GetLiquidResistance, 587
 GetMaxChunkSize_Byte, 587
 GetNumberOfAvailableSamples, 587
 GetPeriod_us, 588
 GetRotaryPositionCode, 588
 GetSampleBufferChunk, 588
 GetSampleRate, 588
 GetSampleVoltageBuffer_uV, 589
 GetScaleFactorU1, 589
 GetScaleFactorU2, 589
 GetUptimeSeconds, 589
 GetWaveform, 589
 IsInternalCalibrationFinished, 590
 IsSamplingFinished, 590
 SetAmplitude_nA, 590
 SetBufferIndex, 590
 SetClampMode, 591
 SetControllerParams, 591
 SetCurrentEnable, 591
 SetExternalLED, 591
 SetLiquidResistance, 592
 SetPeriod_us, 592
 SetWaveform, 592
 StartInternalCalibration, 592
 StartSampling, 592
 StopSampling, 593
CTEERMachineDeviceNet, 593
 ~CTEERMachineDeviceNet, 593
 CTEERMachineDeviceNet, 593
 TEERFunctionNet, 594
CurrentClamp
 Mcs::Usb, 87
CurrentMeasure
 Mcs::Usb, 54
CurrentRangeInNanoAmp

- W2100_StimulusParametersNet, 698
- CurrentResolutionInNanoAmp
 - W2100_StimulusParametersNet, 698
- CUsbDeviceConfigurationFunctionNet, 594
 - !CUsbDeviceConfigurationFunctionNet, 595
 - ~CUsbDeviceConfigurationFunctionNet, 595
 - CUsbDeviceConfigurationFunctionNet, 594
 - SetDeviceId, 595
 - SetDeviceName, 595
- CUsbExceptionNet, 595
 - CUsbExceptionNet, 596
 - Status, 596
- CutoffFrequency
 - CCreateFilterNet, 112
- CVoltageRangeInfoNet
 - CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet, 597
- CW2100_FunctionNet, 597
 - ClearStimulusParametersCache, 599
 - ClearUserDefinedNameCache, 599
 - CW2100_FunctionNet, 599
 - DeselectAllHeadstages, 600
 - DeselectHeadstage, 600
 - GetAccelGyroCurrentRate, 600
 - GetAccelGyroDesiredRate, 600
 - GetAccelGyroEnabled, 600
 - GetAccelRange, 600
 - GetAnalogOutChannel, 600
 - GetAnalogOutFilter, 600
 - GetAudioChannels, 600
 - GetAvailableHeadstages, 600
 - GetBatteryState, 601
 - GetDacRange, 601
 - GetFilterProperties, 601
 - GetFilterProperty, 601
 - GetFPGA FirmwareType, 601
 - GetGyroRange, 601
 - GetHeadstageOnOff, 601
 - GetHeadstageSamplingActive, 601
 - GetMultiHeadstageMode, 601
 - GetPicFirmwareType, 601
 - GetSelectedChannels, 602
 - GetSelectedHeadstageState, 602
 - GetStimulusParametersCache, 602
 - GetStimulusParametersFromSelectedHS, 602
 - GetStimulusParameters, 602
 - GetUserDefinedName, 602
 - GetUserDefinedNameCache, 602
 - GetUserDefinedNameFromSelectedHS, 603
 - PulseGenerator, 604
 - SelectHeadstage, 603
 - SetAccelGyroDesiredRate, 603
 - SetAccelGyroEnabled, 603
 - SetAccelRange, 603
 - SetAnalogOutChannel, 603
 - SetAnalogOutFilter, 603
 - SetAudioChannels, 603
 - SetDacRange, 604
 - SetGyroRange, 604
 - SetHeadstageOnOff, 604
 - SetHeadstageSamplingActive, 604
 - SetHeadstageToSleep, 604
 - SetMultiHeadstageMode, 604
 - SetSelectedChannels, 604
 - Stimulator, 605
- CW2100_FunctionNet::AudioChannelsNet, 92
 - amplification, 92
 - channel, 92
 - dacqgroup, 92
- CW2100_StimulatorFunctionNet, 605
 - BOOST_BIT, 611
 - ClearChannelData, 606
 - CW2100_StimulatorFunctionNet, 606
 - GetBoostAlwaysOnMode, 607
 - GetBoostPreTime, 607
 - GetCurrentRangeInNanoAmp, 607
 - GetCurrentResolutionInNanoAmp, 607
 - GetDACResolution, 607
 - GetDigitalStimulatorTrigger, 607
 - GetDigitalStimulatorTriggerSlope, 608
 - GetNumberOfAnalogChannels, 608
 - GetNumberOfSyncoutChannels, 608
 - GetNumberOfTriggerInputs, 608
 - GetStimulationPatternMemory, 608
 - GetTimeResolutionInNanoSeconds, 608
 - GetTimeSlot, 608
 - GetVoltageRangeInMicroVolt, 608
 - GetVoltageResolutionInMicroVolt, 609
 - GND_SWITCH_BIT, 611
 - PollStatusEvent, 611
 - PrepareData, 609
 - PrepareDataSync, 609
 - SelectTimeSlot, 609
 - SendPreparedData, 609
 - SendStart, 610
 - SendStop, 610
 - SetDigitalStimulatorTrigger, 610
 - SetDigitalStimulatorTriggerSlope, 610
 - StartPoll, 610
 - StopPoll, 610
 - SYNC_BIT0, 611
 - SYNC_BIT1, 611
- CW2100DacqGroupChannelSelectionNet, 611
 - CW2100DacqGroupChannelSelectionNet, 611
- CWarnerUssingDeviceNet, 612
 - !CWarnerUssingDeviceNet, 613
 - ~CWarnerUssingDeviceNet, 612
 - CWarnerUssingDeviceNet, 612
 - WarnerUssingFunction, 613
- CWarnerUssingFunctionNet, 613
 - !CWarnerUssingFunctionNet, 615
 - ~CWarnerUssingFunctionNet, 615
 - CompensateElectrodeOffset, 615
 - CWarnerUssingFunctionNet, 615
 - GetAvailableChambers, 616
 - GetChannelsCountOfChamber, 616

- GetClampMode, 616
- GetComplianceVoltageThreshold, 616
- GetDacPampsPerDigitHighCurrentRange, 617
- GetDacPampsPerDigitLowCurrentRange, 617
- GetDacZero, 617
- GetHighCurrentRange, 618
- GetIdleModeOffset, 618
- GetLiquidResistance, 618
- GetLowCurrentRange, 619
- GetNumberOfAvailableChambers, 619
- GetNumberOfHardwareSlotsForChambers, 619
- GetU1Offset, 619
- GetU1Reference, 621
- GetU2Offset, 621
- GetU2Reference, 621
- GetUnitDescription, 622
- GetUnitExponent, 622
- GetUnitName, 622
- GetUnitsPerDigit, 623
- GetUptimeSeconds, 623
- GetVoltageClampControllerParam_D, 623
- GetVoltageClampControllerParam_I, 624
- GetVoltageClampControllerParam_P, 624
- IsChamberAvailable, 624
- IsHighCurrentMode, 625
- IsInternalCalibrationFinished, 625
- IsPulseEnabled, 625
- SetClampMode, 626
- SetEnablePulse, 626
- SetHighCurrentMode, 626
- SetIdleModeOffset, 627
- SetLiquidResistance, 627
- SetLowCurrentMode, 627
- SetPulse, 627
- SetVoltageClampControllerParam_D, 628
- SetVoltageClampControllerParam_I, 628
- SetVoltageClampControllerParam_P, 628
- WaitForAllChambers, 629
- WaitForChamber, 629
- CWarnerValveControllerDeviceNet, 629
 - !CWarnerValveControllerDeviceNet, 634
 - ~CWarnerValveControllerDeviceNet, 634
- ClearTableName, 634
- ClearValveTable, 634
- CWarnerValveControllerDeviceNet, 634
- GetActiveRunningTableNumber, 634
- GetActiveRunningTableNumberEvent, 651
- GetAnalogThresholdHigh, 634
- GetAnalogThresholdHighEvent, 651
- GetAnalogThresholdLow, 635
- GetAnalogThresholdLowEvent, 651
- GetAnalogVoltage, 635
- GetAnalogVoltageEvent, 651
- GetCurrentEditTableNumber, 635
- GetCurrentNumberOfValves, 635
- GetCurrentNumberOfValvesEvent, 652
- GetDigitalOutPortValve, 636
- GetDigitalOutPortValveEvent, 652
- GetDigitalPortDirection, 636
- GetDigitalPortDirectionEvent, 652
- GetDisplayMode, 636
- GetDisplayModeEvent, 652
- GetTableName, 636
- GetTableNamebyIndex, 637
- GetTableNamebyIndexEvent, 652
- GetTotalNumberOfDigitalPorts, 637
- GetTotalNumberOfTables, 637
- GetTotalNumberOfValves, 637
- GetTotalTableSize, 638
- GetValveActive, 638
- GetValveActiveEvent, 652
- GetValveBoardRevision, 638
- GetValveBoardRevisionEvent, 652
- GetValveBoardRevisionString, 638
- GetValveCurrent, 638
- GetValveDigitalInPort, 639
- GetValveDigitalInPortEvent, 653
- GetValveLedOn, 639
- GetValveLedOnEvent, 653
- GetValveManualGroup, 639
- GetValveManualGroupEvent, 653
- GetValveManualState, 640
- GetValveManualStateEvent, 653
- GetValveMode, 640
- GetValveModeEvent, 653
- GetValvesActiveMap, 640
- GetValvesManualStateMap, 640
- GetValveTableEntry, 641
- IsDigitalOutPortInverted, 641
- IsDigitalOutPortInvertedEvent, 653
- IsValveDigitalInInverted, 641
- IsValveDigitalInInvertedEvent, 653
- IsValveOpen, 642
- IsValveOpenEvent, 654
- IsValveOpenInAnalogMode, 642
- IsValveOpenInAnalogModeEvent, 654
- IsValveOpenInDigitalMode, 642
- IsValveOpenInDigitalModeEvent, 654
- LoadValveTable, 642
- OnGetActiveRunningTableNumber, 643
- OnGetAnalogThresholdHigh, 643
- OnGetAnalogThresholdLow, 643
- OnGetAnalogVoltage, 643
- OnGetCurrentNumberOfValves, 643
- OnGetDigitalOutPortValve, 643
- OnGetDigitalPortDirection, 643
- OnGetDisplayMode, 643
- OnGetTableNamebyIndex, 643
- OnGetValveActive, 644
- OnGetValveBoardRevision, 644
- OnGetValveDigitalInPort, 644
- OnGetValveLedOn, 644
- OnGetValveManualGroup, 644
- OnGetValveManualState, 644
- OnGetValveMode, 644
- OnIsDigitalOutPortInverted, 644

- OnIsValveDigitalInInverted, [645](#)
- OnIsValveOpen, [645](#)
- OnIsValveOpenInAnalogMode, [645](#)
- OnIsValveOpenInDigitalMode, [645](#)
- OnTableEntryChanged, [645](#)
- SetActiveRunningTableNumber, [645](#)
- SetAnalogThresholdHigh, [645](#)
- SetAnalogThresholdLow, [646](#)
- SetCurrentEditTableNumber, [646](#)
- SetDefault, [646](#)
- SetDigitalOutPortInvert, [646](#)
- SetDigitalOutPortValve, [647](#)
- SetDigitalPortDirection, [647](#)
- SetDisplayMode, [647](#)
- SetTableName, [647](#)
- SetTableStep, [648](#)
- SetTableStepAll, [648](#)
- SetValveActive, [648](#)
- SetValveCurrent, [648](#)
- SetValveDigitalInInvert, [649](#)
- SetValveDigitalInPort, [649](#)
- SetValveLedOn, [649](#)
- SetValveManualGroup, [649](#)
- SetValveManualState, [650](#)
- SetValveMode, [650](#)
- SetValvesActiveMap, [650](#)
- SetValvesManualStateMap, [650](#)
- SetValveTableEntry, [650](#)
- StoreValveTable, [651](#)
- TableEntryChangedEvent, [654](#)
- CWarnerValveControllerDeviceTesterFunctionNet, [654](#)
- !CWarnerValveControllerDeviceTesterFunctionNet, [656](#)
- ~CWarnerValveControllerDeviceTesterFunctionNet, [655](#)
- CWarnerValveControllerDeviceTesterFunctionNet, [655](#)
- GetIO, [656](#)
- GetSync, [656](#)
- SetADC, [656](#)
- SetIO, [656](#)
- SetIODirection, [657](#)
- SetTrigger, [657](#)
- SetTriggerSyncDirection, [657](#)
- CWClassicFunctionNet, [657](#)
- CWClassicFunctionNet, [658](#), [659](#)
- GetFilterParametersHeadstage, [659](#)
- GetHasChecksum, [659](#)
- GetHasRedLedHeadstage, [659](#)
- GetHeadstageOnOff, [659](#)
- GetResetFilter, [659](#)
- GetRFConnectionStatus, [659](#)
- GetRFFrequencyHeadstage, [659](#)
- GetRFFrequencyReceiver, [659](#)
- GetRFPower, [660](#)
- GetScanHeadstagesResult, [660](#)
- GetSelectedHeadstage, [660](#)
- GetSerialNumberHeadstage, [660](#)
- GetWPADebugMode, [660](#)
- GetWPAType, [660](#)
- ResetChannelmap, [660](#)
- ScanForHeadstages, [660](#)
- SetChannelmap, [660](#)
- SetFilterParametersHeadstage, [660](#)
- SetHasChecksum, [661](#)
- SetHeadstageOnOff, [661](#)
- SetHWSelectedChannels, [661](#)
- SetResetFilter, [661](#)
- SetRFFrequencyHeadstage, [661](#)
- SetRFFrequencyReceiver, [661](#)
- SetRFFrequencyReceiverEeprom, [661](#)
- SetRFLostBehaviour, [661](#)
- SetRFPower, [662](#)
- SetSelectedHeadstage, [662](#)
- SetSerialNumberHeadstage, [662](#)
- SetWPADebugMode, [662](#)
- SetWPAType, [662](#)
- CWirelessBaseFunctionNet, [662](#)
- CreateWirelessHeadstageSerialNumberString, [663](#)
- CWirelessBaseFunctionNet, [663](#)
- CyclePort
- CMcsUsbNet, [314](#)
- Cypress
- Mcs::Usb, [87](#)
- Cypress_FX1
- Mcs::Usb, [73](#)
- Cypress_FX2
- Mcs::Usb, [73](#)
- Cypress_FX3
- Mcs::Usb, [73](#)
- DAC1Channel
- Mcs::Usb, [56](#)
- DAC2Channel
- Mcs::Usb, [56](#)
- DAC3Channel
- Mcs::Usb, [56](#)
- DAC4Channel
- Mcs::Usb, [56](#)
- DACQ1DigitalGroup
- Mcs::Usb, [57](#)
- DacqDevice
- CSafeISDeviceNet, [491](#)
- dacqgroup
- CW2100_FunctionNet::AudioChannelsNet, [92](#)
- DacqGroupChannelEnumNet
- Mcs::Usb, [57](#)
- DacqMeaGroupTypeEnumNet
- Mcs::Usb, [57](#)
- DacqTrigger
- Mcs::Usb, [62](#)
- DACResolution
- W2100_StimulusParametersNet, [698](#)
- DataModeEnumNet
- Mcs::Usb, [58](#)
- DataState

- HeadStageIDTypeState, [681](#)
- DeepCopy
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [435](#)
- DefineAmplification
 - CPgaDeviceNet, [400](#)
- DefineFrequencyRange
 - CPgaDeviceNet, [400](#)
- DefineNumAmplifications
 - CPgaDeviceNet, [400](#)
- DefineNumFrequencyRanges
 - CPgaDeviceNet, [400](#)
- DeselectAllHeadstages
 - CW2100_FunctionNet, [600](#)
- DeselectHeadstage
 - CW2100_FunctionNet, [600](#)
- DEST_FX3_TARGET_MASK
 - Mcs::Usb, [54](#)
- DEST_TARGET1
 - Mcs::Usb, [53](#)
- DEST_TARGET10
 - Mcs::Usb, [53](#)
- DEST_TARGET11
 - Mcs::Usb, [53](#)
- DEST_TARGET12
 - Mcs::Usb, [53](#)
- DEST_TARGET13
 - Mcs::Usb, [54](#)
- DEST_TARGET14
 - Mcs::Usb, [54](#)
- DEST_TARGET15
 - Mcs::Usb, [54](#)
- DEST_TARGET2
 - Mcs::Usb, [53](#)
- DEST_TARGET3
 - Mcs::Usb, [53](#)
- DEST_TARGET4
 - Mcs::Usb, [53](#)
- DEST_TARGET5
 - Mcs::Usb, [53](#)
- DEST_TARGET6
 - Mcs::Usb, [53](#)
- DEST_TARGET7
 - Mcs::Usb, [53](#)
- DEST_TARGET8
 - Mcs::Usb, [53](#)
- DEST_TARGET9
 - Mcs::Usb, [53](#)
- DEST_TARGET_MASK
 - Mcs::Usb, [54](#)
- DetectChipType
 - CCMOSMea_FunctionNet, [99](#)
- DEVICE_NOT_FOUND
 - Mcs::Usb, [63](#)
- DeviceArrival
 - CMcsUsbListNet, [307](#)
- DeviceData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- DeviceDataLength
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- DeviceEnumNet
 - Mcs::Usb, [58](#)
- DeviceHasNoHeadstage
 - Mcs::Usb, [65](#), [78](#)
- DeviceId
 - CMcsUsbListEntryNet, [304](#)
- DeviceIdNet, [663](#)
 - BcdDevice, [664](#)
 - BusType, [664](#)
 - DeviceIdNet, [664](#)
 - IdProduct, [664](#)
 - IdVendor, [664](#)
 - operator=, [664](#)
- DeviceName
 - CMcsUsbListEntryNet, [304](#)
- DeviceNotConnected
 - Mcs::Usb, [64](#), [65](#), [78](#)
- DeviceRemoval
 - CMcsUsbListNet, [307](#)
- DeviceRunStatus
 - Mcs::Usb, [61](#), [70](#), [79](#), [83](#), [89](#)
- DigDataFromReceiver
 - Mcs::Usb, [89](#)
- Digital
 - Mcs::Usb, [73](#), [90](#)
- DigitalData
 - Mcs::Usb, [61](#), [70](#), [79](#), [83](#), [89](#)
- DigitalDatastreamEnableEnumNet
 - Mcs::Usb, [59](#)
- DigitalGroup
 - Mcs::Usb, [58](#)
- DigitalIn
 - Mcs::Usb, [59](#), [60](#), [70](#), [79](#), [82](#), [89](#)
- DigitalInOfOutPort
 - Mcs::Usb, [60](#), [70](#), [79](#), [82](#), [89](#)
- DigitalInPort
 - Mcs::Usb, [55](#)
- DigitalInReserverd
 - Mcs::Usb, [59](#)
- DigitalMux
 - Mcs::Usb, [55](#)
- DigitalOut
 - Mcs::Usb, [59](#)
- DigitalOutReg
 - Mcs::Usb, [55](#)
- DigitalOutReserved
 - Mcs::Usb, [60](#)
- DigitalOutStimulator
 - Mcs::Usb, [61](#), [70](#), [79](#), [83](#), [89](#)
- DigitalPulse
 - Mcs::Usb, [60](#), [70](#), [79](#), [82](#), [89](#)
- DigitalReg
 - Mcs::Usb, [56](#)

- DigitalSource
 - DigitalSource< digitalsourceenum >, 665
- DigitalSource< digitalsourceenum >, 665
 - DigitalSource, 665
 - MaxBitNumber, 665
 - MaxBitNumberStatic, 665
 - size, 665
 - Source, 666
- DigitalSourceEnumNet
 - Mcs::Usb, 60
- DigitalSourceGeneral, 666
 - DigitalSourceGeneral, 666
 - MaxBitNumber, 666, 667
 - size, 667
 - Source, 667
- DigitalStimulatorTriggerEventEnumNet
 - Mcs::Usb, 61
- DigitalStimulatorTriggerSlopeEnumNet
 - Mcs::Usb, 61
- DigitalTargetEnumNet
 - Mcs::Usb, 61
- Digout
 - Mcs::Usb, 62
- DigOutStim
 - Mcs::Usb, 60
- DigOutStimulatorStartTrigger
 - Mcs::Usb, 62
- DigOutStimulatorStopTrigger
 - Mcs::Usb, 62
- Digstream
 - Mcs::Usb, 62
- DigStreamFromReceiver
 - Mcs::Usb, 89
- DigStreamToReceiver
 - Mcs::Usb, 62
- Dilutor
 - Mcs::Usb, 76
- DisableAutoReset
 - CStg200xDownloadBasicNet, 544
- DisableMultiFileMode
 - CStg200xDownloadNet, 552
- Disconnect
 - CMcsUsbNet, 314
- DisConnectDevice
 - CRadioControlledDevicesNet, 434
- DongleS
 - Mcs::Usb, 74
- DoRamp
 - CRoboDacqNet, 449
- Dotriapot
 - Mcs::Usb, 74
- DoubleToInt
 - Mcs::Usb, 63
- DownloadFirmware
 - CMcsUsbFactoryNet, 292
- DownloadOnly
 - Mcs::Usb, 81
- DriverVersionNet, 667
 - ~DriverVersionNet, 668
 - DriverVersionNet, 668
 - DriverVersionNet::FormatVersion, 668
 - GetDestinationCode, 669
 - GetDestinationName, 670
 - GetMajor, 670
 - GetMinor, 671
 - GetNumEntries, 671
 - GetSerialNumber, 671
 - GetStatus, 672
 - GetVersionInt, 672
 - GetVersionString, 673
- DriverVersionNet::FormatVersion
 - DriverVersionNet, 668
- DSP
 - FirmwareDestinationNames, 674
 - Mcs::Usb, 51
- DSPAnalogGroup
 - Mcs::Usb, 69, 79
- DSPDataGroup
 - Mcs::Usb, 57, 88
- DSPDigitalGroup
 - Mcs::Usb, 69, 79
- DummyCommand
 - CLIH3DeviceNet, 190
- Duration
 - CStimulusFunctionNet::SidebandData, 694
- eCube
 - Mcs::Usb, 75, 78
- eCubeHeadstage
 - Mcs::Usb, 65
- ElectricalStimulation
 - HeadStageIDType, 677
- ElectrodeDacMuxEnumNet
 - Mcs::Usb, 62
- ElectrodeModeEnumNet
 - Mcs::Usb, 62
- ElectrodeOffset
 - Mcs::Usb, 87
- emAutomatic
 - Mcs::Usb, 62
- emManual
 - Mcs::Usb, 62
- EmptyKey
 - CMcsUsbNet, 314
- Emu_GetCellCapacity
 - CRoboDacqNet, 450
- Emu_GetCellPotential
 - CRoboDacqNet, 450
- Emu_GetCellResists
 - CRoboDacqNet, 450
- Emu_GetElectrodeResists
 - CRoboDacqNet, 450
- Emu_GetNoise
 - CRoboDacqNet, 450
- Emu_SetCellCapacity
 - CRoboDacqNet, 450
- Emu_SetCellPotential

- CRoboDacqNet, [450](#)
- Emu_SetCellResists
 - CRoboDacqNet, [450](#)
- Emu_SetElectrodeResists
 - CRoboDacqNet, [450](#)
- Emu_SetNoise
 - CRoboDacqNet, [450](#)
- EnableAnalogOut
 - CSCUFunctionNet, [495](#)
- EnableAutoReset
 - CStg200xDownloadBasicNet, [544](#)
- EnableChannelsInGroup
 - CCMOSMea_FunctionNet, [99](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, [116](#)
- EnableChecksum
 - CMeaDeviceNet, [345](#)
 - COctoPotDeviceNet, [389](#)
- EnableDigitalIn
 - CMeaDeviceNet, [346](#)
 - COctoPotDeviceNet, [389](#)
- EnableExceptions
 - CMcsUsbNet, [314](#)
- EnableMultiFileMode
 - CStg200xDownloadNet, [552](#)
- EnableQueue
 - CRoboDeviceNet, [465](#)
- EnableTimestamp
 - CMeaDeviceNet, [347](#)
 - COctoPotDeviceNet, [389](#)
- EnableUserTrigger
 - CLIH3DeviceNet, [190](#)
- Encapsulator
 - Mcs::Usb, [76](#)
- enCMosMeaChipType
 - Mcs::Usb, [62](#)
- EnSTG200x_STATUS
 - Mcs::Usb, [63](#)
- Entry
 - HeadStageIDType, [678](#)
- EOAndCRC
 - Mcs::Usb, [56](#)
- Equals
 - CMcsUsbListEntryNet, [301](#)
 - HeadStageIDType, [678](#)
 - HeadstageIDTypeObject, [680](#)
- EraseEepromRegisterPreconfig
 - CMcsUsbNet, [314](#)
- EraseFilterParameterPermanent
 - CFilterConfigurationNet, [132](#)
 - CFilterConfigurationRegisterNet, [133](#), [134](#)
- ErasePermanentAdcOffset
 - CLIH3DeviceNet, [191](#)
- ErasePermanentDacOffset
 - CLIH3DeviceNet, [191](#)
- Error_Callback_Aquisition_Stopped
 - CMcsUsbDacqNet, [286](#)
- Error_Callback_Data_Lost
 - CMcsUsbDacqNet, [286](#)
- Error_Callback_Frames_Lost
 - CMcsUsbDacqNet, [286](#)
- Error_Callback_Packet_Error
 - CMcsUsbDacqNet, [286](#)
- Error_Callback_Queue_Full
 - CMcsUsbDacqNet, [287](#)
- Error_Callback_RingQueue_Full
 - CMcsUsbDacqNet, [287](#)
- ErrorEvent
 - CMcsUsbDacqNet, [287](#)
- ExternBCTester
 - Mcs::Usb, [74](#)
- ExternDTester
 - Mcs::Usb, [74](#)
- ExternSTester
 - Mcs::Usb, [74](#)
- FactoryReset
 - CTcxDeviceNet, [572](#)
- Falling
 - Mcs::Usb, [61](#)
- FCB
 - Mcs::Usb, [74](#)
- FCX
 - Mcs::Usb, [74](#)
- Feedback
 - Mcs::Usb, [60](#), [70](#), [79](#), [82](#), [89](#)
- FeedbackGetSampleTimerCount
 - CMeaFeedbackFunctionNet, [354](#)
- FeedbackHigh
 - Mcs::Usb, [60](#)
- FeedbackLow
 - Mcs::Usb, [60](#)
- FeedbackReg
 - Mcs::Usb, [56](#)
- FeedbackSetAnalogSource
 - CMeaFeedbackFunctionNet, [354](#)
- FeedbackSetChannelFilter
 - CMeaFeedbackFunctionNet, [354](#)
- FeedbackSetDigitalMapping
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetFeedback
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetFilterOff
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetFilterParameter
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetFilterParameter32
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetGlobalChannelFilter
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetIIRFilterParameter
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetLogic
 - CMeaFeedbackFunctionNet, [355](#)
- FeedbackSetMkFilter

- CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetNumberOfLogics
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetNumberOfRateCounter
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetNumberOfRateDetectors
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetNumberOfSpikeDetectors
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetNumberOfTriggers
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetRateCounter
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetRateDetector
 - CMeaFeedbackFunctionNet, [356](#)
- FeedbackSetSpikeDetectorThreshold
 - CMeaFeedbackFunctionNet, [357](#)
- FeedbackSetTrigger
 - CMeaFeedbackFunctionNet, [357](#)
- FilterActive
 - CFilterPropertyNet, [136](#)
- FilterAttributeEnumNet
 - Mcs::Usb, [63](#)
- FilterBand
 - CFilterPropertyNet, [136](#)
- FilterBandEnumNet
 - Mcs::Usb, [63](#)
- FilterCalculationDirectionEnumNet
 - Mcs::Usb, [63](#)
- FilterFamily
 - CFilterPropertyNet, [136](#)
- FilterFamilyEnumNet
 - Mcs::Usb, [63](#)
- FilterType
 - CFilterPropertyNet, [136](#)
- FilterTypeEnumNet
 - Mcs::Usb, [64](#)
- FindEndpoints
 - CGenericDevelopDeviceNet, [152](#)
- FindFilter
 - CCreateFilterNet, [111](#)
- FindFirmwareVersionMagicInBuffer
 - CMcsUsbFactoryNet, [292](#)
- FindReference
 - CRoboDeviceNet, [465](#)
- FindReferencel
 - CRoboStatorDeviceNet, [484](#)
- FindReferencePhase0
 - CRoboDeviceNet::RoboMainLowLevelCommands, [686](#)
- FindReferencePhase0XY
 - CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands, [692](#)
- FindReferenceXY
 - CRoboStatorDeviceNet, [484](#)
- FindReferenceZ
 - CRoboStatorDeviceNet, [484](#)
- Finished
 - Mcs::Usb, [81](#)
- FirePressurePulse
 - CPPCFunctionNet, [416](#)
- FirmwareDestinationNames, [673](#)
 - Altera, [674](#)
 - Bootstrap, [674](#)
 - BUS1_MCSBUS1, [674](#)
 - BUS1_MCSBUS2, [674](#)
 - DSP, [674](#)
 - FPGA2, [674](#)
 - FPGA3, [674](#)
 - FPGA4, [674](#)
 - FPGA5, [674](#)
 - FPGA6, [674](#)
 - MCSBUS1, [675](#)
 - MCSBUS10, [675](#)
 - MCSBUS11, [675](#)
 - MCSBUS12, [675](#)
 - MCSBUS13, [675](#)
 - MCSBUS2, [675](#)
 - MCSBUS3, [675](#)
 - MCSBUS4, [675](#)
 - MCSBUS5, [675](#)
 - MCSBUS6, [675](#)
 - MCSBUS7, [675](#)
 - MCSBUS8, [676](#)
 - MCSBUS9, [676](#)
 - MCU1, [676](#)
 - PIC, [676](#)
 - PIC2, [676](#)
 - PIC3, [676](#)
 - PIC4, [676](#)
 - USB, [676](#)
- FluidControlDevice
 - CSafeISDeviceNet, [491](#)
- ForceStatusEvent
 - CStg200xDownloadBasicNet, [544](#)
 - CStimulusFunctionNet, [559](#)
- FPGA10
 - Mcs::Usb, [52](#)
- FPGA10_BASE
 - Mcs::Usb, [53](#)
- FPGA10_GOLD
 - Mcs::Usb, [53](#)
- FPGA11
 - Mcs::Usb, [52](#)
- FPGA11_BASE
 - Mcs::Usb, [53](#)
- FPGA11_GOLD
 - Mcs::Usb, [53](#)
- FPGA12
 - Mcs::Usb, [52](#)
- FPGA12_BASE
 - Mcs::Usb, [53](#)
- FPGA12_GOLD
 - Mcs::Usb, [53](#)
- FPGA13
 - Mcs::Usb, [52](#)

FPGA13_BASE
Mcs::Usb, [53](#)
FPGA13_GOLD
Mcs::Usb, [53](#)
FPGA14
Mcs::Usb, [52](#)
FPGA14_BASE
Mcs::Usb, [53](#)
FPGA14_GOLD
Mcs::Usb, [53](#)
FPGA15
Mcs::Usb, [52](#)
FPGA15_BASE
Mcs::Usb, [53](#)
FPGA15_GOLD
Mcs::Usb, [53](#)
FPGA16
Mcs::Usb, [53](#)
FPGA16_BASE
Mcs::Usb, [53](#)
FPGA16_GOLD
Mcs::Usb, [53](#)
FPGA2
FirmwareDestinationNames, [674](#)
Mcs::Usb, [52](#)
FPGA2_BASE
Mcs::Usb, [53](#)
FPGA2_GOLD
Mcs::Usb, [53](#)
FPGA3
FirmwareDestinationNames, [674](#)
Mcs::Usb, [52](#)
FPGA3_BASE
Mcs::Usb, [53](#)
FPGA3_GOLD
Mcs::Usb, [53](#)
FPGA4
FirmwareDestinationNames, [674](#)
Mcs::Usb, [52](#)
FPGA4_BASE
Mcs::Usb, [53](#)
FPGA4_GOLD
Mcs::Usb, [53](#)
FPGA5
FirmwareDestinationNames, [674](#)
Mcs::Usb, [52](#)
FPGA5_BASE
Mcs::Usb, [53](#)
FPGA5_GOLD
Mcs::Usb, [53](#)
FPGA6
FirmwareDestinationNames, [674](#)
Mcs::Usb, [52](#)
FPGA6_BASE
Mcs::Usb, [53](#)
FPGA6_GOLD
Mcs::Usb, [53](#)
FPGA7
Mcs::Usb, [52](#)
FPGA7_BASE
Mcs::Usb, [53](#)
FPGA7_GOLD
Mcs::Usb, [53](#)
FPGA8
Mcs::Usb, [52](#)
FPGA8_BASE
Mcs::Usb, [53](#)
FPGA8_GOLD
Mcs::Usb, [53](#)
FPGA9
Mcs::Usb, [52](#)
FPGA9_BASE
Mcs::Usb, [53](#)
FPGA9_GOLD
Mcs::Usb, [53](#)
FPGA_BASE
Mcs::Usb, [53](#)
FPGA_BOOTSTRAP
Mcs::Usb, [53](#)
FPGA_GOLD
Mcs::Usb, [53](#)
FPGA_NORMAL
Mcs::Usb, [51](#)
FpgalEnumNet
Mcs::Usb, [64](#)
FrameContextGroup
Mcs::Usb, [58](#)
FromIntPtr
StgStatusNet, [695](#)
FromPtr
StgStatusNet, [695](#)
FullCharge
Mcs::Usb, [67](#)
FullSpeed
Mcs::Usb, [69](#)
FunkDongleS
Mcs::Usb, [74](#)
FX3MCSDDataAddress
CMcsUsbFactoryNet, [297](#)
FX3MCSDDataDeviceIdOffset
CMcsUsbFactoryNet, [297](#)
FX3MCSDDataIFB1ImageOffset
CMcsUsbFactoryNet, [297](#)
FX3MCSDDataIFB2ImageOffset
CMcsUsbFactoryNet, [297](#)
FX3MCSDDataVersionOffset
CMcsUsbFactoryNet, [297](#)
FYIProgram
CFYIDeviceNet, [144](#)
FYITemp
CFYIDeviceNet, [144](#)
Gain
CMeaDeviceNet, [350](#)
Gated_High_Active
Mcs::Usb, [76](#)
Gated_Low_Active

- Mcs::Usb, [76](#)
- GE2100
 - Mcs::Usb, [75](#)
- Generic
 - Mcs::Usb, [74](#)
- Get2AnalogInput
 - CMcsBus_SensorNet, [223](#)
- Get2DigitalInput
 - CMcsBus_SensorNet, [223](#)
- Get4ADC
 - CMcsBus_SensorNet, [223](#)
- Get4ADCAverage
 - CMcsBus_SensorNet, [223](#)
- Get4ADCCatchampAverageShift
 - CMcsBus_SensorNet, [223](#)
- Get4ADCMode
 - CMcsBus_SensorNet, [223](#)
- Get4DAC
 - CMcsBus_SensorNet, [223](#)
- GetAbsMaxCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [383](#)
- GetAccelGyroCurrentRate
 - CW2100_FunctionNet, [600](#)
- GetAccelGyroDesiredRate
 - CW2100_FunctionNet, [600](#)
- GetAccelGyroEnabled
 - CW2100_FunctionNet, [600](#)
- GetAccelRange
 - CW2100_FunctionNet, [600](#)
- GetActiveRunningTableNumber
 - CWarnerValveControllerDeviceNet, [634](#)
- GetActiveRunningTableNumberEvent
 - CWarnerValveControllerDeviceNet, [651](#)
- GetAdapterCode
 - CMealImpedanceDeviceNet, [358](#)
 - CTEERFunctionNet, [585](#)
- GetAdapterType
 - CMcsUsbDacqNet, [262](#)
- GetAdc
 - CFluidControlDeviceNet, [138](#)
- GetAdcDataFormat
 - CMcsUsbDacqNet, [262](#)
- GetADCInputOffset
 - CCMOSMea_FunctionNet, [99](#)
- GetAdcOffset
 - CLIH3DeviceNet, [191](#)
 - COctoPotDeviceNet, [389](#)
- GetAdcOffsetU1
 - CTEERFunctionNet, [585](#)
- GetAdcOffsetU2
 - CTEERFunctionNet, [585](#)
- GetADCs
 - CMcsBus_SensorNet, [223](#)
- GetADCsLoop
 - CMcsBus_SensorNet, [224](#)
- GetAdcZero
 - CMcsUsbDacqNet, [263](#)
- GetAirpressure
 - CRoboDeviceNet, [465](#)
- GetAirpressureLimit
 - CRoboDeviceNet, [465](#)
- GetAirValve
 - CRoboDeviceNet, [465](#)
- GetAllDigout
 - CRoboDacqNet, [451](#)
- GetAmplification
 - CPgaDeviceNet, [400](#)
- GetAmplitude_nA
 - CTEERFunctionNet, [585](#)
- GetAnalogGain
 - CMeaDeviceNet, [347](#)
- GetAnalogOutADCRange
 - CSCUFunctionNet, [495](#)
- GetAnalogOutChannel
 - CW2100_FunctionNet, [600](#)
- GetAnalogOutChannels
 - CSCUFunctionNet, [495](#)
- GetAnalogOutDACRange
 - CSCUFunctionNet, [496](#)
- GetAnalogOutFilter
 - CW2100_FunctionNet, [600](#)
- GetAnalogRanges
 - CStg200xBasicNet, [513](#)
- GetAnalogResolution
 - CStg200xBasicNet, [513](#)
- GetAnalogThresholdHigh
 - CWarnerValveControllerDeviceNet, [634](#)
- GetAnalogThresholdHighEvent
 - CWarnerValveControllerDeviceNet, [651](#)
- GetAnalogThresholdLow
 - CWarnerValveControllerDeviceNet, [635](#)
- GetAnalogThresholdLowEvent
 - CWarnerValveControllerDeviceNet, [651](#)
- GetAnalogValueUnit
 - CMcsUsbDacqNet, [263](#)
- GetAnalogVoltage
 - CPPCFunctionNet, [416](#)
 - CPPS_FunctionNet, [425](#)
 - CWarnerValveControllerDeviceNet, [635](#)
- GetAnalogVoltageEvent
 - CWarnerValveControllerDeviceNet, [651](#)
- GetAnalogVoltageRange
 - CPPCFunctionNet, [416](#)
- GetAnalogVoltages
 - CPPS_FunctionNet, [425](#)
- GetArraySize
 - CMealImpedanceDeviceNet, [358](#)
- GetAudioChannels
 - CMeaAudioFunctionNet, [332](#)
 - CW2100_FunctionNet, [600](#)
- GetAudioOutDacParameter
 - CLIH3DeviceNet, [191](#)
- GetAutocalibrationDisabled
 - CStg200xBasicNet, [513](#)
- GetAvailableBaseSamplerates
 - CCMOSMeaDeviceNet, [108](#)

- CGrapheneASICDeviceNet, 165
- GetAvailableChambers
 - CWarnerUssingFunctionNet, 616
- GetAvailableDeviceList
 - CRFFunctionNet, 440
- GetAvailableDeviceListEx
 - CRFFunctionNet, 440
- GetAvailableHeadstages
 - CSCUFunctionNet, 496
 - CW2100_FunctionNet, 600
- GetAvailableHeadstagesEvent
 - CSCUFunctionNet, 506
- GetAvailableMemory
 - CStg200xBasicNet, 515
 - CStimulusFunctionNet, 560
- GetAvailableSampleRates
 - CMcsUsbDacqNet::CHWInfo, 180
- GetAvailableStateList
 - CRFFunctionNet, 440
- GetAvailableStateListEx
 - CRFFunctionNet, 440
- GetAvailableVoltageRangesInMicroVolt
 - CMcsUsbDacqNet::CHWInfo, 180
- GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt
 - CMcsUsbDacqNet::CHWInfo, 181
- GetAxisConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, 686
- GetAxisLED
 - CRoboocyte2DeviceNet, 482
- GetAxisParametersSignedEeprom
 - CMcsBus_AxisParametersNet, 198
- GetAxisParametersUnsignedEeprom
 - CMcsBus_AxisParametersNet, 198
- GetBaseFrequency
 - CRFFunctionNet, 441
- GetBaseSamplerate
 - CCMOSMeaDeviceNet, 109
- GetBath
 - CCMOSMea_FunctionNet, 99
- GetBathMode
 - CCMOSMea_FunctionNet, 99
- GetBatteryState
 - CW2100_FunctionNet, 601
- GetBatteryVoltage
 - CMultiBatteryChargerDeviceNet, 366
- GetBiQuad
 - CCreateFilterNet, 112
- GetBiQuads
 - CCreateFilterNet, 112
- GetBlankingEnable
 - CStg200xBasicNet, 515
- GetBoardTemp
 - CTcxDeviceNet, 572
- GetBoostAlwaysOnMode
 - CW2100_StimulatorFunctionNet, 607
- GetBoostPreTime
 - CW2100_StimulatorFunctionNet, 607
- GetBubbleState
 - CPPS_FunctionNet, 425
- GetBubbleStatus
 - CMcsBus_SensorNet, 224
- GetBuffer
 - CGenericDevelopDeviceNet, 152
- GetBusAddress
 - CMcsBusNet, 237
- GetBusAddressEeprom
 - CMcsBusNet, 237
- GetByteBuffer
 - CGenericDevelopDeviceNet, 153
- GetBytesAvailable
 - CSerialPortNet, 507
- GetBytesPerSample
 - CTEERFunctionNet, 586
- GetCalibration
 - CTcxDeviceNet, 572
- GetCalibrationDecp
 - CTcxDeviceNet, 572
- GetCalibrationMax
 - CTcxDeviceNet, 572
- GetCalibrationMin
 - CTcxDeviceNet, 572
- GetCapacityC
 - CRoboDacqNet, 451
- GetCapacityV
 - CRoboDacqNet, 451
- GetCapacityX
 - CRoboDacqNet, 451
- GetCardinalDacqSamplerate
 - CInterfaceboardFunctionNet, 187
- GetCardinalStgOutputrate
 - CInterfaceboardFunctionNet, 187
- GetChannel
 - CSw2to64DeviceNet, 568
- GetChannelDataFillSize
 - CMcsUsbDacqNet, 263
- GetChannelDataI16
 - CCMOSMeaDeviceNet, 109
- GetChannelDataI32
 - CCMOSMeaDeviceNet, 109
- GetChannelDataUI16
 - CCMOSMeaDeviceNet, 109
- GetChannelDataUI32
 - CCMOSMeaDeviceNet, 109
- GetChannelLayout
 - CMcsUsbDacqNet, 263
- GetChannels
 - CMultiBatteryChargerDeviceNet, 366
 - CSw2to64DeviceNet, 568
- GetChannelsCountOfChamber
 - CWarnerUssingFunctionNet, 616
- GetChannelsInBlock
 - CMcsUsbDacqNet, 263
- GetChannelState
 - CMultiBatteryChargerDeviceNet, 366
- GetChargeCapacity

- CMultiBatteryChargerDeviceNet, 367
- GetChargeCurrent
 - CMultiBatteryChargerDeviceNet, 367
- GetChargingMode
 - CMultiBatteryChargerDeviceNet, 367
- GetChargingPCoefficient
 - CMultiBatteryChargerDeviceNet, 368
- GetChecksumFromFX3Image
 - CMcsUsbFactoryNet, 292
- GetCheckVoltage
 - CokuvisionStimulatorDeviceNet, 392
- GetClampAmpSerialNumber
 - CRoboDacqNet, 451
- GetClampMode
 - CTEERFunctionNet, 586
 - CWarnerUssingFunctionNet, 616
- GetCMOSDataDictionary
 - CCMOSMeaDeviceNet, 109
- GetCoilCommunication
 - CPositionIIDeviceNet, 403
- GetColorRgb
 - CMultiwellOptoStimFunctionNet, 383
- GetColorStr
 - CMultiwellOptoStimFunctionNet, 383
- GetCommand
 - CMcsBusNet, 237, 238
 - CPedoterDeviceNet, 397
 - CRoboDacqNet, 451
- GetComplianceVoltageThreshold
 - CWarnerUssingFunctionNet, 616
- GetConfiguration
 - CMcsUsbNet, 315
- GetConfigurationBit
 - CRoboDacqNet, 451
- GetConfigurationBitAxc
 - CRoboDacqNet, 451
- GetConfigurationBitBlu_Led
 - CRoboDacqNet, 451
- GetConfigurationBitBlu_LedToggleFast
 - CRoboDacqNet, 451
- GetConfigurationBitBlu_LedToggleSlow
 - CRoboDacqNet, 451
- GetConfigurationBitCC_Gen
 - CRoboDacqNet, 452
- GetConfigurationBitCV_Gen
 - CRoboDacqNet, 452
- GetConfigurationBitRC_Gen
 - CRoboDacqNet, 452
- GetConfigurationBitRed_Led
 - CRoboDacqNet, 452
- GetConfigurationBitRed_LedSaturation
 - CRoboDacqNet, 452
- GetConfigurationBitRed_LedToggleFast
 - CRoboDacqNet, 452
- GetConfigurationBitRed_LedToggleSlow
 - CRoboDacqNet, 452
- GetConfigurationBitRelais
 - CRoboDacqNet, 452
- GetConfigurationBitRV_Gen
 - CRoboDacqNet, 452
- GetConfigurationBits
 - CRoboDacqNet, 452
- GetConfigurationBitStream
 - CRoboDacqNet, 452
- GetConfigurationBitSupply
 - CRoboDacqNet, 453
- GetConnectedDevice
 - CRFFFunctionNet, 441
- GetControllerParams
 - CTEERFunctionNet, 586
- GetCrossTalkOffset
 - CRoboDacqNet, 453
- GetCrossTalkOptimum
 - CRoboDacqNet, 453
- GetCur2VolResistance
 - CGrapheneFunctionNet, 168
- GetCurrent
 - CTcxDeviceNet, 572
- GetCurrentAirvalve
 - CRoboDeviceNet, 465
- GetCurrentAirvalveLimit
 - CRoboDeviceNet, 465
- GetCurrentCycle
 - CMeaCoatDeviceNet, 339
- GetCurrentEditTableNumber
 - CWarnerValveControllerDeviceNet, 635
- GetCurrentEnable
 - CTEERFunctionNet, 586
- GetCurrentFactor
 - CokuvisionStimulatorDeviceNet, 392
- GetCurrentNumberOfValves
 - CWarnerValveControllerDeviceNet, 635
- GetCurrentNumberOfValvesEvent
 - CWarnerValveControllerDeviceNet, 652
- GetCurrentPosition
 - CRoboDeviceNet, 466
- GetCurrentPositionI
 - CRoboStatorDeviceNet, 484
- GetCurrentPositionXY
 - CRoboStatorDeviceNet, 484
- GetCurrentPositionZ
 - CRoboStatorDeviceNet, 485
- GetCurrentRangeByIndex
 - CStg200xBasicNet, 516
- GetCurrentRangeInNanoAmp
 - CStg200xBasicNet, 516
 - CStimulusFunctionNet, 560
 - CW2100_StimulatorFunctionNet, 607
- GetCurrentResolutionInNanoAmp
 - CStg200xBasicNet, 516
 - CStimulusFunctionNet, 560
 - CW2100_StimulatorFunctionNet, 607
- GetCycle
 - CMeaCleanDeviceNet, 335
- GetCycles
 - CMeaCleanDeviceNet, 335

- CMeaCoatDeviceNet, 339
- GetD
 - CTcxDeviceNet, 573
- GetDacAmplificationFactor
 - CStg200xBasicNet, 517
- GetDacIdleValue
 - CLIH3DeviceNet, 192
- GetDACOffset
 - CGrapheneFunctionNet, 169
 - COkuvisionStimulatorDeviceNet, 393
- GetDacOffset
 - CDacCalibrationFunctionNet, 114
 - CLIH3DeviceNet, 192
 - COctoPotDeviceNet, 389
- GetDacPampsPerDigitHighCurrentRange
 - CWarnerUssingFunctionNet, 617
- GetDacPampsPerDigitLowCurrentRange
 - CWarnerUssingFunctionNet, 617
- GetDacqRunStatus
 - CLIH3DeviceNet, 192
- GetDacRange
 - CW2100_FunctionNet, 601
- GetDACResolution
 - CStg200xBasicNet, 517
 - CStimulusFunctionNet, 560
 - CW2100_StimulatorFunctionNet, 607
- GetDACs
 - CMcsBus_SensorNet, 224
- GetDacUseldleValue
 - CLIH3DeviceNet, 192
- GetDacZero
 - CTEERFunctionNet, 587
 - CWarnerUssingFunctionNet, 617
- GetDataFormat
 - CMcsUsbDacqNet, 263
- GetDataMode
 - CMcsUsbDacqNet, 263
- GetDDecp
 - CTcxDeviceNet, 573
- GetDebugData
 - CPositionIIDeviceNet, 403
- GetDestination
 - CMcsUsbFactoryNet, 292
- GetDestinationCode
 - DriverVersionNet, 669
- GetDestinationDisplayLabel
 - CMcsUsbFactoryNet, 292
- GetDestinationName
 - CMcsUsbFactoryNet, 292, 293
 - DriverVersionNet, 670
- GetDestinationSerialNumber
 - CMcsUsbFactoryNet, 293
- GetDestinationTargetAddress
 - CMcsUsbFactoryNet, 293
- GetDetectionThreshold
 - CMcsBus_SensorNet, 224
- GetDetectorValue
 - CMcsBus_SensorNet, 224
- GetDevice
 - CTcxDeviceNet, 573
- GetDeviceCannotStallOutRequests
 - CMcsUsbNet, 315
- GetDeviceCapableSpeed
 - CMcsUsbNet, 315
- GetDeviceEnum
 - CMcsUsbNet, 315
- GetDeviceGroupChannelInfos
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, 116, 117
- GetDeviceId
 - CMcsUsbNet, 315
- GetDeviceList
 - CPositionImpDeviceNet, 411
- GetDeviceNames
 - CRadioControlledDevicesNet, 434
- GetDeviceRootHubVendorEnum
 - CMcsUsbNet, 315
- GetDeviceRootHubVendorID
 - CMcsUsbNet, 315
- GetDeviceRootHubVendorName
 - CMcsUsbNet, 315
- GetDeviceSpeed
 - CMcsUsbNet, 316
- GetDeviceType
 - CTcxDeviceNet, 573
- GetDevname
 - CTcxDeviceNet, 573
- GetDigin
 - CFluidControlDeviceNet, 138
- GetDiginState
 - CLIH3DeviceNet, 193
- GetDiginValue
 - CStg200xBasicNet, 517
- GetDigitalData
 - CMeaDigitalDataFunctionNet, 352
- GetDigitalIn
 - CPPCFunctionNet, 418
 - CPPS_FunctionNet, 425
- GetDigitalOutPortValve
 - CWarnerValveControllerDeviceNet, 636
- GetDigitalOutPortValveEvent
 - CWarnerValveControllerDeviceNet, 652
- GetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, 636
- GetDigitalPortDirectionEvent
 - CWarnerValveControllerDeviceNet, 652
- GetDigitalSource
 - CMcsUsbDacqNet, 264, 265
- GetDigitalStimulatorTrigger
 - CW2100_StimulatorFunctionNet, 607
- GetDigitalStimulatorTriggerSlope
 - CW2100_StimulatorFunctionNet, 608
- GetDigout
 - CFluidControlDeviceNet, 139

- CRoboDacqNet, 453
- GetDigoutMode
 - CStg200xBasicNet, 517
- GetDigoutValue
 - CStg200xBasicNet, 517
- GetDIO
 - CMcsBus_FYIExtensionNet, 201
- GetDischargeCapacity
 - CMultiBatteryChargerDeviceNet, 368
- GetDischargeCurrent
 - CMultiBatteryChargerDeviceNet, 368
- GetDischargeCurrentSetPoint
 - CMultiBatteryChargerDeviceNet, 368
- GetDisplayMode
 - CWarnerValveControllerDeviceNet, 636
- GetDisplayModeEvent
 - CWarnerValveControllerDeviceNet, 652
- GetDisplayText
 - CRoboDacqNet, 453
- GetDMax
 - CTcxDeviceNet, 573
- GetDMin
 - CTcxDeviceNet, 573
- GetDownsampleFactor
 - CRoboDacqNet, 453
- GetDSPHighPassByIndex
 - CIntanMea_FunctionNet, 183
- GetDuration
 - CMeaCoatDeviceNet, 340
- GetDuty
 - CTcxDeviceNet, 573
- GetEEPromPage
 - CLIH3DeviceNet, 193
- GetElectrodeDacMux
 - CStg200xBasicNet, 518
- GetElectrodeEnable
 - CStg200xBasicNet, 518, 519
- GetElectrodeMode
 - CStg200xBasicNet, 519, 520
- GetEnableAmplifierProtectionSwitch
 - CStg200xBasicNet, 520
- GetEnabledChannelsInGroup
 - CCMOSMea_FunctionNet, 100
 - CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, 117
- GetEnableHeaterLimit
 - CTcxDeviceNet, 574
- GetEnableThermocouple
 - CTcxDeviceNet, 574
- GetEntry
 - CMcsUsbListEntryNet, 302
- GetEntryCount
 - CMcsUsbListEntryNet, 303
- GetEnumerationSpeed
 - CMeaDeviceNet, 347
- GetErrorAirpressure
 - CRoboDeviceNet, 466
- GetErrorCurrentAirvalve
 - CRoboDeviceNet, 466
- GetErrorMessage
 - CMcsUsbDacqNet, 266
- GetErrorText
 - CMcsUsbNet, 316
- GetErrorVoltage12V
 - CRoboDeviceNet, 466
- GetErrorVoltage5V
 - CRoboDeviceNet, 466
- GetErrorVoltageAirvalve
 - CRoboDeviceNet, 466
- GetErrorVoltageRs485A
 - CRoboDeviceNet, 466
- GetErrorVoltageRs485B
 - CRoboDeviceNet, 466
- GetErrorVoltageValves
 - CRoboDeviceNet, 467
- GetEventData
 - CPositionIIDeviceNet, 404
- GetExternalElectrodeEnable
 - CStg200xBasicNet, 521
- GetFAAmplification
 - CStg200xBasicNet, 521
- GetFilter
 - CRoboDacqNet, 453
- GetFilterAttributes
 - CFilterConfigurationNet, 132
- GetFilterCoeffs
 - CRoboDacqNet, 453
- GetFilterParametersHeadstage
 - CWClassicFunctionNet, 659
- GetFilterProperties
 - CSCUFunctionNet, 496
 - CW2100_FunctionNet, 601
- GetFilterProperty
 - CMcsUsbDacqNet, 266
 - CSCUFunctionNet, 497
 - CW2100_FunctionNet, 601
- GetFinalDischargeVoltage
 - CMultiBatteryChargerDeviceNet, 369
- GetFirmwareVersion
 - CMcsUsbNet, 316
- GetFirmwareVersionFromFile
 - CMcsUsbFactoryNet, 293
- GetFirmwareVersionFromHexFile
 - CMcsUsbFactoryNet, 293
- GetFPGA FirmwareType
 - CW2100_FunctionNet, 601
- GetFrameErrorCounter
 - CTEERFunctionNet, 587
- GetFrequency
 - CRadioControlledDevicesNet, 435
- GetFrequencyRange
 - CPgaDeviceNet, 401
- GetGain
 - CMeaDeviceNet, 347

- CPgaDeviceNet, [401](#)
- GetGate
 - CCMOSMea_FunctionNet, [100](#)
- GetGilsonDevice
 - CRobococyte2DeviceNet, [482](#)
- GetGlobalRepeat
 - CDigOutStimulatorFunctionNet, [123](#)
- GetGNDI
 - CCMOSMea_FunctionNet, [100](#)
- GetGroupADCBits
 - CCMOSMea_FunctionNet, [100](#)
- GetGroupChannelBitmaskBySelect
 - CCMOSMea_FunctionNet, [100](#)
- GetGroupChannelBitmaskHS1NCBathCurrent
 - CCMOSMea_FunctionNet, [100](#), [101](#)
- GetGroupChannelBitmaskHS1NCCol2Current
 - CCMOSMea_FunctionNet, [101](#)
- GetGroupChannelBitmaskHS1NChipTemp
 - CCMOSMea_FunctionNet, [101](#)
- GetGroupChannelBitmaskHS1Sidebands
 - CCMOSMea_FunctionNet, [101](#)
- GetGroupChannelBitmaskHS1TriggerStatus
 - CCMOSMea_FunctionNet, [101](#), [102](#)
- GetGroupChannelBitmaskIFDigChannels
 - CCMOSMea_FunctionNet, [102](#)
- GetGroupChannelBitmaskInterfaceADC
 - CCMOSMea_FunctionNet, [102](#)
- GetGroupChannelBitmaskPacketFrameContext
 - CCMOSMea_FunctionNet, [102](#)
- GetGroupChannelBitmaskSTG1DACSignal
 - CCMOSMea_FunctionNet, [102](#), [103](#)
- GetGroupChannelData16
 - CMcsUsbDacqNet, [266](#)
- GetGroupChannelData32
 - CMcsUsbDacqNet, [266](#)
- GetGroupChannelDataUI16
 - CMcsUsbDacqNet, [267](#)
- GetGroupChannelDataUI32
 - CMcsUsbDacqNet, [267](#)
- GetGroupDCOffset
 - CCMOSMea_FunctionNet, [103](#)
- GetGroupID
 - CCMOSMea_FunctionNet, [103](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-qGroupChannelEnumTemplateNet, Dac-qGroupChannelEnumTemplate, CDevice-GroupChannelInfoTemplateNet >, [117](#)
- GetGroupNumberOfChannels
 - CCMOSMea_FunctionNet, [103](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-qGroupChannelEnumTemplateNet, Dac-qGroupChannelEnumTemplate, CDevice-GroupChannelInfoTemplateNet >, [117](#)
- GetGroupResolutionPerDigit
 - CCMOSMea_FunctionNet, [103](#)
- GetGroupSampleSize
 - CCMOSMea_FunctionNet, [104](#)
- CDacqGroupChannelSelectionTemplateNet< Dac-qGroupChannelEnumTemplateNet, Dac-qGroupChannelEnumTemplate, CDevice-GroupChannelInfoTemplateNet >, [117](#)
- GetGroupType
 - CCMOSMea_FunctionNet, [104](#)
 - CDacqGroupChannelSelectionTemplateNet< Dac-qGroupChannelEnumTemplateNet, Dac-qGroupChannelEnumTemplate, CDevice-GroupChannelInfoTemplateNet >, [118](#)
- GetGroupUnit
 - CCMOSMea_FunctionNet, [104](#)
- GetGyroRange
 - CW2100_FunctionNet, [601](#)
- GetHardwareMaxRange
 - CMcsUsbDacqNet, [268](#)
- GetHardwareMinRange
 - CMcsUsbDacqNet, [268](#)
- GetHardwareRevision
 - CMcsUsbNet, [316](#)
- GetHasChecksum
 - CWClassicFunctionNet, [659](#)
- GetHashCode
 - HeadstageIDTypeObject, [680](#)
- GetHasRedLedHeadstage
 - CWClassicFunctionNet, [659](#)
- GetHasThermocouple
 - CTcxDeviceNet, [574](#)
- GetHeadstage
 - CStg200xBasicNet, [521](#)
- GetHeadstageActive
 - CMcsUsbNet, [317](#)
- GetHeadstageAdcBits
 - CSCUFunctionNet, [497](#)
- GetHeadstageAdcRangeInMicroVolt
 - CSCUFunctionNet, [497](#)
- GetHeadstageDacBits
 - CSCUFunctionNet, [498](#)
- GetHeadstageDacCurrentRangeInMicroAmpere
 - CSCUFunctionNet, [498](#)
- GetHeadstageDacCurrentResolutionInNanoAmpere
 - CSCUFunctionNet, [498](#)
- GetHeadstageDacVoltageRangeInMilliVolt
 - CSCUFunctionNet, [499](#)
- GetHeadstageDacVoltageResolutionInMicroVolt
 - CSCUFunctionNet, [499](#)
- GetHeadstageGainInPer mille
 - CSCUFunctionNet, [499](#)
- GetHeadstageID
 - CMcsUsbNet, [317](#)
 - CSCUFunctionNet, [500](#)
- GetHeadstageNumberOfAnalogChannels
 - CSCUFunctionNet, [500](#)
- GetHeadstageNumberOfStimulationChannels
 - CSCUFunctionNet, [500](#)
- GetHeadstageOnOff
 - CW2100_FunctionNet, [601](#)
 - CWClassicFunctionNet, [659](#)

- GetHeadstagePowerStateAtStart
 - CSCUFunctionNet, [501](#)
- GetHeadstagePresent
 - CMcsUsbNet, [317](#)
- GetHeadstageSamplerate
 - CSCUFunctionNet, [501](#)
- GetHeadstageSamplingActive
 - CW2100_FunctionNet, [601](#)
- GetHeadstageSerialNumber
 - CSCUFunctionNet, [501](#)
- GetHeaterLimit
 - CTcxDeviceNet, [574](#)
- GetHeaterTemp
 - CTcxDeviceNet, [574](#)
- GetHighCurrentRange
 - CWarnerUssingFunctionNet, [618](#)
- GetHighpassFilterEnable
 - CFilterConfigurationNet, [132](#)
- GetHWConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, [686](#)
- GetHWRevision
 - CRoboDeviceNet::RoboMainLowLevelCommands, [686](#)
- GetHWRevisionEeprom
 - CMcsBusNet, [238](#)
- GetI
 - CTcxDeviceNet, [574](#)
- GetIC
 - CRoboDacqNet, [453](#)
- GetIClamp
 - CRoboDacqNet, [453](#)
- GetICoeff
 - CRobo_FYITemp_FunctionNet, [445](#)
- GetICOffset
 - CRoboDacqNet, [453](#)
- GetIDecp
 - CTcxDeviceNet, [574](#)
- GetIdleModeOffset
 - CWarnerUssingFunctionNet, [618](#)
- GetIGain
 - CRoboDacqNet, [454](#)
- GetIMax
 - CTcxDeviceNet, [575](#)
- GetIMin
 - CTcxDeviceNet, [575](#)
- GetImpedanceResult
 - CIntanMea_FunctionNet, [183](#)
- GetImpedanceTestFrequency
 - CMealImpedanceDeviceNet, [358](#)
- GetImpId
 - CPositionImpDeviceNet, [412](#)
- GetImplantCurrentSetpoint
 - CPositionIIDeviceNet, [404](#)
- GetImplantResult
 - CPositionIIDeviceNet, [405](#)
- GetImplantState
 - CPositionIIDeviceNet, [405](#)
- GetInMovement
 - CRoboDeviceNet, [467](#)
- GetIntanRegister
 - CIntanMea_FunctionNet, [183](#)
- GetIntBuffer
 - CGenericDevelopDeviceNet, [153](#)
- GetIO
 - CWarnerValveControllerDeviceTesterFunctionNet, [656](#)
- GetIOOut
 - CTcxDeviceNet, [575](#)
- GetIOVoltage
 - CInterfaceboard2FunctionNet, [185](#)
- GetLastAnswer
 - CGilsonDeviceNet, [164](#)
- GetLastUSBError
 - CMcsUsbNet, [317](#)
- GetLatency
 - CMcsBus_SensorNet, [224](#)
- GetLatencyCounter
 - CMcsBus_SensorNet, [224](#)
- GetLayoutConfiguration
 - CMEA2100x256FunctionNet, [330](#)
- GetLEDSwitch
 - CMcsBus_ExtensionNet, [200](#)
- GetLength
 - CRobo_FYIProgram_FunctionNet, [443](#)
- GetLiquidResistance
 - CTEERFunctionNet, [587](#)
 - CWarnerUssingFunctionNet, [618](#)
- GetListmodeIndexRange
 - CStg200xBasicNet, [521](#)
- GetListmodeTriggerSource
 - CStg200xBasicNet, [522](#)
- GetLowCurrentRange
 - CWarnerUssingFunctionNet, [619](#)
- GetLowerFrequencyByIndex
 - CIntanMea_FunctionNet, [183](#)
- GetMajor
 - DriverVersionNet, [670](#)
- GetMaxChunkSize_Byte
 - CTEERFunctionNet, [587](#)
- GetMaxCurrent
 - CMeaCoatDeviceNet, [340](#)
- GetMaxDurationHighCurrentInMicroSec
 - CMultiwellOptoStimFunctionNet, [384](#)
- GetMaxDutyCycleHighCurrent
 - CMultiwellOptoStimFunctionNet, [384](#)
- GetMaxHeaterPowerMultiwell
 - CTcxDeviceNet, [575](#)
- GetMaxNoPressure
 - CRoboDeviceNet::RoboMainLowLevelCommands, [686](#)
- GetMaxNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [686](#)
- GetMaxNumberOfHeadstages
 - CSCUFunctionNet, [502](#)

GetMaxNumOfColumns
 CCMOSMea_FunctionNet, 104

GetMaxP
 CTcxDeviceNet, 575

GetMaxpDecp
 CTcxDeviceNet, 575

GetMaxpMax
 CTcxDeviceNet, 575

GetMaxpMin
 CTcxDeviceNet, 576

GetMaxPower
 CokuvisionStimulatorDeviceNet, 393
 CRobo_FYITemp_FunctionNet, 445

GetMaxPressureWaitTime
 CRoboDeviceNet::RoboMainLowLevelCommands, 686

GetMaxReadableColumns
 CCMOSMeaDeviceNet, 109

GetMaxSamplingFrequency
 CMcsUsbDacqNet, 268

GetMaxStimulusChannelsPerHeadstage
 CSCUFunctionNet, 502

GetMaxVoltage
 CMeaCleanDeviceNet, 335
 CokuvisionStimulatorDeviceNet, 393

GetMCAcceleration
 CMcsBus_MotorControlNet, 205

GetMCAccelerationEeprom
 CMcsBus_MotorControlNet, 205

GetMCAccelerationShortCommand
 CMcsBus_MotorControlNet, 205

GetMCAxisRevisionEeprom
 CMcsBus_MotorControlNet, 205

GetMCBreakCurrent
 CMcsBus_MotorControlNet, 205

GetMCBreakCurrentEeprom
 CMcsBus_MotorControlNet, 206

GetMCConfig
 CMcsBus_MotorControlNet, 206

GetMCConfigEeprom
 CMcsBus_MotorControlNet, 206

GetMCCurrent
 CMcsBus_MotorControlNet, 206

GetMCCurrentEeprom
 CMcsBus_MotorControlNet, 206

GetMCCurrentMode
 CMcsBus_MotorControlNet, 206

GetMCCurrentModeEeprom
 CMcsBus_MotorControlNet, 206

GetMCCurrentModeShortCommand
 CMcsBus_MotorControlNet, 207

GetMCCurrentPosition
 CMcsBus_MotorControlNet, 207

GetMCCurrentShortCommand
 CMcsBus_MotorControlNet, 207

GetMCCurrentSpeed
 CMcsBus_MotorControlNet, 207

GetMCMaxAcceleration
 CMcsBus_MotorControlNet, 207

GetMCMaxAccelerationEeprom
 CMcsBus_MotorControlNet, 207

GetMCMaxCurrent
 CMcsBus_MotorControlNet, 207

GetMCMaxCurrentEeprom
 CMcsBus_MotorControlNet, 208

GetMCMaxSpeed
 CMcsBus_MotorControlNet, 208

GetMCMaxSpeedEeprom
 CMcsBus_MotorControlNet, 208

GetMCMaxTravel
 CMcsBus_MotorControlNet, 208

GetMCMaxTravelEeprom
 CMcsBus_MotorControlNet, 208

GetMCMaxTravelShortCommand
 CMcsBus_MotorControlNet, 208

GetMCMovement
 CMcsBus_MotorControlNet, 208

GetMCNewPosition
 CMcsBus_MotorControlNet, 209

GetMCOutputOnOff
 CMcsBus_MotorControlNet, 209

GetMCPhase
 CMcsBus_MotorControlNet, 209

GetMCPhaseOffset
 CMcsBus_MotorControlNet, 209

GetMCReference
 CMcsBus_MotorControlNet, 209

GetMCReferenceCurrent
 CMcsBus_MotorControlNet, 209

GetMCReferenceCurrentEeprom
 CMcsBus_MotorControlNet, 209

GetMCRegulatorGain
 CMcsBus_MotorControlNet, 210

GetMCRegulatorGainEeprom
 CMcsBus_MotorControlNet, 210

GetMcsBus_Extension
 CRoboocyte2DeviceNet, 482

GetMCScalingFactor
 CMcsBus_MotorControlNet, 210

GetMCScalingFactorEeprom
 CMcsBus_MotorControlNet, 210

GetMCSpeed
 CMcsBus_MotorControlNet, 210

GetMCSpeedEeprom
 CMcsBus_MotorControlNet, 210

GetMCSpeedShortCommand
 CMcsBus_MotorControlNet, 210

GetMCSpeedUnitEeprom
 CMcsBus_MotorControlNet, 211

GetMCStandbyCurrent
 CMcsBus_MotorControlNet, 211

GetMCStandbyCurrentEeprom
 CMcsBus_MotorControlNet, 211

GetMCStandbyTime
 CMcsBus_MotorControlNet, 211

GetMCStandbyTimeEeprom

- CMcsBus_MotorControlNet, 211
- GetMea21UsbPort
 - CMcsUsbNet, 318
- GetMeaLayout
 - CMcsUsbDacqNet, 268
- GetMemoryUsageDAC
 - CStg200xDownloadBasicNet, 544
- GetMemoryUsageSyncout
 - CStg200xDownloadBasicNet, 544
- GetMinimalThreshold
 - CMcsBus_SensorNet, 224
- GetMinNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, 686
- GetMinor
 - DriverVersionNet, 671
- GetMinPressure
 - CRoboDeviceNet, 467
 - CRoboDeviceNet::RoboMainLowLevelCommands, 686
- GetMinPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetMinSamplingFrequencyStepsize
 - CMcsUsbDacqNet, 269
- GetMinVoltage
 - CMeaCleanDeviceNet, 336
- GetModeSelect
 - CPulseGeneratorFunctionNet, 432
- GetModuleCurrent
 - CStg200xDownloadNet, 552
- GetModuleTemp
 - CStg200xDownloadNet, 552
- GetMovementError
 - CRoboDeviceNet, 467
- GetMovePump
 - CMcsBus_SensorNet, 225
- GetMultiHeadstageMode
 - CW2100_FunctionNet, 601
- GetMultiplexedDataChannelsInBlock
 - CStimulusFunctionNet, 561
- GetNanoVoltsPerKelvin
 - CMcsBus_TempSensorNet, 230
- GetNeurochipMemoryData
 - CCMOSMea_FunctionNet, 104
- GetNeurochipMemorySize
 - CCMOSMea_FunctionNet, 105
- GetNIC_MS
 - CRoboDacqNet, 454
- GetNUC_MS
 - CRoboDacqNet, 454
- GetNumAmplifications
 - CPgaDeviceNet, 401
- GetNumber
 - CMeaSwitchDeviceNet, 361
 - CSw2to64DeviceNet, 568
- GetNumberOfAnalogChannels
 - CStg200xBasicNet, 522
- CStimulusFunctionNet, 561
- CW2100_StimulatorFunctionNet, 608
- GetNumberOfAudioChannels
 - CMeaAudioFunctionNet, 332
- GetNumberOfAvailableChambers
 - CWarnerUssingFunctionNet, 619
- GetNumberOfAvailableSamples
 - CTEERFunctionNet, 587
- GetNumberOfChannels
 - CDigOutStimulatorFunctionNet, 124
- GetNumberOfDataBits
 - CMcsUsbDacqNet, 269
- GetNumberOfDevices
 - CMcsUsbListNet, 306
- GetNumberOfHardwareSlotsForChambers
 - CWarnerUssingFunctionNet, 619
- GetNumberOfHWADCCChannels
 - CMcsUsbDacqNet::CHWInfo, 181
- GetNumberOfHWDACPaths
 - CStg200xBasicNet, 522
- GetNumberOfHWDigitalChannels
 - CMcsUsbDacqNet::CHWInfo, 181
- GetNumberOfStimulationElectrodes
 - CStg200xBasicNet, 523
- GetNumberOfStimulationSourcesPerElectrode
 - CStg200xBasicNet, 523
- GetNumberOfSupportedGroups
 - CCMOSMea_FunctionNet, 105
- CDacqGroupChannelSelectionTemplateNet< Dac-
qGroupChannelEnumTemplateNet, Dac-
qGroupChannelEnumTemplate, CDevice-
GroupChannelInfoTemplateNet >, 118
- GetNumberOfSyncoutChannels
 - CStg200xBasicNet, 523
 - CW2100_StimulatorFunctionNet, 608
- GetNumberOfTriggerInputs
 - CStg200xBasicNet, 523
 - CW2100_StimulatorFunctionNet, 608
- GetNumConfigurations
 - CMcsUsbNet, 318
- GetNumControlChannels
 - CTcxDeviceNet, 576
- GetNumDestinations
 - CMcsUsbFactoryNet, 293
- GetNumDevices
 - CTcxDeviceNet, 576
- GetNumEntries
 - DriverVersionNet, 671
- GetNumFrequencyRanges
 - CPgaDeviceNet, 401
- GetNumMeasureChannels
 - CTcxDeviceNet, 576
- GetNUV_MS
 - CRoboDacqNet, 454
- GetOffsetCurrent
 - CMeaCoatDeviceNet, 340
- GetOnOff
 - CPositionIIDeviceNet, 405

- CTcxDeviceNet, 576
- GetOutputCurrent
 - CMeaCoatDeviceNet, 340
- GetOutputRate
 - CStg200xBasicNet, 523
- GetOutputVoltage
 - CMeaCleanDeviceNet, 336
- GetP
 - CTcxDeviceNet, 576
- GetParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetPattern
 - CMeaSwitchDeviceNet, 361
- GetPatternBool
 - CMeaSwitchDeviceNet, 361
- GetPauseDuration
 - CMeaCoatDeviceNet, 340
- GetPCoeff
 - CRobo_FYITemp_FunctionNet, 445
- GetPDecp
 - CTcxDeviceNet, 576
- GetPeriod
 - CPulseGeneratorFunctionNet, 432
- GetPeriod_us
 - CTEERFunctionNet, 588
- GetPermanentCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, 384
- GetPGain
 - CRoboDacqNet, 454
- GetPhases
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetPicFirmwareType
 - CW2100_FunctionNet, 601
- GetPiezoState
 - CMcsBus_SensorNet, 225
- GetPlateClampLockState
 - CMultiwellDeviceNet, 376
- GetPlateClampState
 - CMultiwellDeviceNet, 376
- GetPlateClampStateByHeadstage
 - CMultiwellCallbackFunctionNet, 373
- GetPlateClampStateByHeadstageEvent
 - CMultiwellCallbackFunctionNet, 374
- GetPlateMux
 - CMultiwellDeviceNet, 377
- GetPlateType
 - CMultiwellDeviceNet, 377
- GetPMax
 - CTcxDeviceNet, 576
- GetPMin
 - CTcxDeviceNet, 576
- GetPoti
 - CMcsUsbDacqNet, 269
- GetPOut
 - CTcxDeviceNet, 576
- GetPowerMuxPlate
 - CMultiwellDeviceNet, 378
- GetPowerStrength
 - CPositionIIDeviceNet, 406
- GetPressure
 - CMcsBus_SensorNet, 225
- GetPressureOffset
 - CMcsBus_SensorNet, 225
- GetPressureRange
 - CPPCFunctionNet, 418
- GetPulseform
 - COkuvisionStimulatorDeviceNet, 393
- GetPulseLength
 - CPulseGeneratorFunctionNet, 432
- GetPumpCouple
 - CPPS_FunctionNet, 425
- GetPumpEnableSpeedRatio
 - CPPS_FunctionNet, 425
- GetPumpFastOnOff
 - CPPS_FunctionNet, 425
- GetPumpFastSpeed
 - CPPS_FunctionNet, 425
- GetPumpFunctionSpeeds
 - CPPS_FunctionNet, 426
- GetPumpManualOnOff
 - CPPS_FunctionNet, 426
- GetPumpMaxSpeed
 - CPPS_FunctionNet, 426
- GetPumpModeType
 - CPPCFunctionNet, 418
 - CPPS_FunctionNet, 426
- GetPumpSpeed
 - CRoboFluidDeviceNet, 478
- GetPumpSpeedRatio
 - CPPS_FunctionNet, 426
- GetPumpSpeedUnit
 - CPPCFunctionNet, 419
 - CPPS_FunctionNet, 426
- GetPWM
 - CFluidControlDeviceNet, 139
- GetPwrOut
 - CTcxDeviceNet, 577
- GetPwrSet
 - CTcxDeviceNet, 577
- GetRatedCapacity
 - CMultiBatteryChargerDeviceNet, 369
- GetReady
 - CMealImpedanceDeviceNet, 358
- GetRecordingNumber
 - CRoboDacqNet, 454
- GetReferenceElectrodeMode
 - CSCUFunctionNet, 502
- GetReferenceElectrodeSwitchState
 - CSCUFunctionNet, 502
- GetReferenceTemperature
 - CFluidControlDeviceNet, 139
- GetRegulationTimeouts
 - CMcsBus_SensorNet, 225
- GetRegulatorFactor

- CMcsBus_SensorNet, 226
- GetRegulatorOnOff
 - CMcsBus_SensorNet, 226
 - CRobo_FYITemp_FunctionNet, 445
- GetRegulatorStatus
 - CMcsBus_SensorNet, 226
- GetRepeats
 - CProgramPressureCurveNet, 430
- GetRes1
 - CTcxDeviceNet, 577
- GetRes2
 - CTcxDeviceNet, 577
- GetResetFilter
 - CWClassicFunctionNet, 659
- GetResistanceC
 - CRoboDacqNet, 454
- GetResistanceV
 - CRoboDacqNet, 454
- GetResolutionPerDigit
 - CMcsUsbDacqNet, 269
- GetResS
 - CTcxDeviceNet, 577
- GetResult
 - CMealImpedanceDeviceNet, 358
- GetResX
 - CTcxDeviceNet, 577
- GetRFConnectionStatus
 - CWClassicFunctionNet, 659
- GetRFFrequency
 - CPositionImpDeviceNet, 412
- GetRFFrequencyHeadstage
 - CWClassicFunctionNet, 659
- GetRFFrequencyReceiver
 - CWClassicFunctionNet, 659
- GetRFPower
 - CWClassicFunctionNet, 660
- GetRoboDacq
 - CRoboocyte2DeviceNet, 482
- GetRoboFluidDevice
 - CEncapsulatorDeviceNet, 127
 - CRoboocyte2DeviceNet, 482
- GetRotaryPositionCode
 - CTEERFunctionNet, 588
- GetRotatePump
 - CMcsBus_SensorNet, 226
- GetROut
 - CTcxDeviceNet, 577
- GetRTC
 - CokuvisionStimulatorDeviceNet, 393
 - CPositionIIDeviceNet, 406
- GetSampleBufferChunk
 - CTEERFunctionNet, 588
- GetSampleInterval
 - CLIH3DeviceNet, 193
- GetSamplePeriode
 - CMcsBus_SensorNet, 226
- GetSampleRate
 - CTEERFunctionNet, 588
- GetSamplerate
 - CMcsUsbDacqNet, 269
- GetSampleVoltageBuffer_uV
 - CTEERFunctionNet, 589
- GetScaleFactorU1
 - CTEERFunctionNet, 589
- GetScaleFactorU2
 - CTEERFunctionNet, 589
- GetScanHeadstagesResult
 - CWClassicFunctionNet, 660
- GetScreen
 - CRoboDacqNet, 454
- GetSearchReferenceFastAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetSearchReferenceFastSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetSearchReferenceFineAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetSearchReferenceFineSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, 687
- GetSearchReferenceMethod
 - CRoboDeviceNet::RoboMainLowLevelCommands, 688
- GetSearchReferenceMoveOut
 - CRoboDeviceNet::RoboMainLowLevelCommands, 688
- GetSearchReferenceOffsetPos
 - CRoboDeviceNet::RoboMainLowLevelCommands, 688
- GetSelectedChannels
 - CW2100_FunctionNet, 602
- GetSelectedHeadstage
 - CWClassicFunctionNet, 660
- GetSelectedHeadstageState
 - CW2100_FunctionNet, 602
- GetSensorType
 - CTcxDeviceNet, 578
- GetSerialNumber
 - CMcsUsbNet, 318
 - DriverVersionNet, 671
- GetSerialNumberHeadstage
 - CWClassicFunctionNet, 660
- GetSetpoint
 - CTcxDeviceNet, 578
- GetSetpointDecp
 - CTcxDeviceNet, 578
- GetSetpointMax
 - CTcxDeviceNet, 578
- GetSetpointMin
 - CTcxDeviceNet, 578
- GetShortBuffer
 - CGenericDevelopDeviceNet, 154
- GetSimulation
 - CRoboDacqNet, 454

GetSingleHeater
 CMcsBus_FYIExtensionNet, 201

GetSingleValve
 CFluidControlDeviceNet, 139
 CRoboFluidDeviceNet, 478

GetSlope
 CMeaCleanDeviceNet, 336
 CMeaCoatDeviceNet, 341

GetSoftwareKey
 CMcsUsbNet, 318

GetSoftwareKeyString
 CMcsUsbNet, 318

GetSollPressure
 CMcsBus_SensorNet, 226

GetSollTemp
 CRobo_FYITemp_FunctionNet, 445

GetSourceBulk
 CCMOSMea_FunctionNet, 105

GetSourceDrain
 CCMOSMea_FunctionNet, 105

GetSourceGate
 CCMOSMea_FunctionNet, 105

GetStartTriggerSlope
 CDigOutStimulatorFunctionNet, 124

GetState
 CRFFFunctionNet, 441
 CRobo_FYIProgram_FunctionNet, 444

GetStateDebugData
 CPositionIIDeviceNet, 406

GetStateEventData
 CPositionIIDeviceNet, 407

GetStatus
 CMcsUsbNet, 318
 DriverVersionNet, 672

GetStatusOfLastCommand
 CMcsUsbNet, 318

GetStgProgramInfo
 CStg200xBasicNet, 523, 524

GetStgVersionInfo
 CStg200xBasicNet, 524

GetStimulationPatternMemory
 CW2100_StimulatorFunctionNet, 608

GetStimulatorStatus
 CokuvisionStimulatorDeviceNet, 393

GetStimulusParametersCache
 CW2100_FunctionNet, 602

GetStimulusParametersFromSelectedHS
 CW2100_FunctionNet, 602

GetStimulusSites
 CCMOSMea_FunctionNet, 105

GetStimulusParameters
 CW2100_FunctionNet, 602

GetStopTriggerSlope
 CDigOutStimulatorFunctionNet, 124

GetSubChannel
 CMcsBus_MotorControlNet, 211

GetSupplyVoltage
 CPPCFunctionNet, 419
 CPPS_FunctionNet, 426

GetSweepCount
 CStg200xDownloadBasicNet, 545

GetSync
 CWarnerValveControllerDeviceTesterFunctionNet, 656

GetSyncoutMap
 CStg200xBasicNet, 524

GetSyncState
 CMcsBus_SensorNet, 226

GetTableName
 CWarnerValveControllerDeviceNet, 636

GetTableNamebyIndex
 CWarnerValveControllerDeviceNet, 637

GetTableNamebyIndexEvent
 CWarnerValveControllerDeviceNet, 652

GetTablepointer
 CRetinaLedDeviceNet, 437

GetTemperatur
 CMcsBus_TempSensorNet, 230

GetTestMode
 CRFFFunctionNet, 441

GetThermocoupleCalibration
 CFluidControlDeviceNet, 140
 CTcxDeviceNet, 578

GetThermocoupleNanovoltPerKelvin
 CFluidControlDeviceNet, 140
 CTcxDeviceNet, 578

GetThermocoupleReferenceTemp
 CTcxDeviceNet, 579

GetThermocoupleTemp
 CTcxDeviceNet, 579

GetThermocoupleTempAbs
 CTcxDeviceNet, 579

GetThermocoupleTemperature
 CFluidControlDeviceNet, 140

GetThermoOffset
 CMcsBus_TempSensorNet, 230

GetThermoTemp
 CMcsBus_TempSensorNet, 231

GetThermoVoltage
 CMcsBus_TempSensorNet, 231

GetTimeInPause
 CMeaCoatDeviceNet, 341

GetTimeInPlateau
 CMeaCoatDeviceNet, 341

GetTimeResolutionInNanoSeconds
 CW2100_StimulatorFunctionNet, 608

GetTimeSlot
 CW2100_StimulatorFunctionNet, 608

GetTotalMemory
 CStg200xBasicNet, 525
 CStimulusFunctionNet, 561

GetTotalNumberOfDigitalPorts
 CWarnerValveControllerDeviceNet, 637

GetTotalNumberOfTables
 CWarnerValveControllerDeviceNet, 637

GetTotalNumberOfValves

- CWarnerValveControllerDeviceNet, 637
- GetTotalTableSize
 - CWarnerValveControllerDeviceNet, 638
- GetTouchPadEnable
 - CMultiwellDeviceNet, 378
- GetTrigger
 - CStg200xDownloadBasicNet, 545
- GetTriggerSource
 - CStg200xBasicNet, 525
- GetU1Offset
 - CWarnerUssingFunctionNet, 619
- GetU1Reference
 - CWarnerUssingFunctionNet, 621
- GetU2Offset
 - CWarnerUssingFunctionNet, 621
- GetU2Reference
 - CWarnerUssingFunctionNet, 621
- GetUByteBuffer
 - CGenericDevelopDeviceNet, 155
- GetUC
 - CRoboDacqNet, 454
- GetUClamp
 - CRoboDacqNet, 455
- GetUCOffset
 - CRoboDacqNet, 455
- GetUIntA
 - CFilterCoefficientsNet, 130
- GetUIntB
 - CFilterCoefficientsNet, 130
- GetUIntBuffer
 - CGenericDevelopDeviceNet, 155
- GetUnit
 - CTcxDeviceNet, 579
- GetUnitDescription
 - CWarnerUssingFunctionNet, 622
- GetUnitExponent
 - CWarnerUssingFunctionNet, 622
- GetUnitName
 - CWarnerUssingFunctionNet, 622
- GetUnitsPerDigit
 - CWarnerUssingFunctionNet, 623
- GetUOut
 - CTcxDeviceNet, 579
- GetUpdateDisplay
 - CRoboDacqNet, 455
- GetUpperFrequencyByIndex
 - CIntanMea_FunctionNet, 184
- GetUptimeSeconds
 - CTEERFunctionNet, 589
 - CWarnerUssingFunctionNet, 623
- GetUSBDeviceIDFromFX3Image
 - CMcsUsbFactoryNet, 294
- GetUsbListEntries
 - CMcsUsbListNet, 306
- GetUsbListEntry
 - CMcsUsbListNet, 306
 - CMcsUsbNet, 319
- GetUseBubble
 - CPPS_FunctionNet, 426
- GetUserCodeFromBitFile
 - CMcsUsbFactoryNet, 294
- GetUserCodeFromFlash
 - CMcsUsbFactoryNet, 294
- GetUserDefinedName
 - CW2100_FunctionNet, 602
- GetUserDefinedNameCache
 - CW2100_FunctionNet, 602
- GetUserDefinedNameFromSelectedHS
 - CW2100_FunctionNet, 603
- GetUserParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, 688
- GetUShortBuffer
 - CGenericDevelopDeviceNet, 156
- GetUV
 - CRoboDacqNet, 455
- GetUVOffset
 - CRoboDacqNet, 455
- GetValue
 - CTcxDeviceNet, 579
- GetValueHires
 - CTcxDeviceNet, 580
- GetValve
 - CFluidControlDeviceNet, 141
 - CRoboFluidDeviceNet, 478
- GetValve1
 - CRobo_FYIProgram_FunctionNet, 444
- GetValve2
 - CRobo_FYIProgram_FunctionNet, 444
- GetValveActive
 - CPPCFunctionNet, 419
 - CWarnerValveControllerDeviceNet, 638
- GetValveActiveEvent
 - CWarnerValveControllerDeviceNet, 652
- GetValveBoardRevision
 - CWarnerValveControllerDeviceNet, 638
- GetValveBoardRevisionEvent
 - CWarnerValveControllerDeviceNet, 652
- GetValveBoardRevisionString
 - CWarnerValveControllerDeviceNet, 638
- GetValveCurrent
 - CWarnerValveControllerDeviceNet, 638
- GetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, 639
- GetValveDigitalInPortEvent
 - CWarnerValveControllerDeviceNet, 653
- GetValveLedOn
 - CWarnerValveControllerDeviceNet, 639
- GetValveLedOnEvent
 - CWarnerValveControllerDeviceNet, 653
- GetValveManualGroup
 - CWarnerValveControllerDeviceNet, 639
- GetValveManualGroupEvent
 - CWarnerValveControllerDeviceNet, 653
- GetValveManualState
 - CWarnerValveControllerDeviceNet, 640

GetValveManualStateEvent
 CWarnerValveControllerDeviceNet, 653

GetValveMode
 CWarnerValveControllerDeviceNet, 640

GetValveModeEvent
 CWarnerValveControllerDeviceNet, 653

GetValves
 CMcsBus_FYIExtensionNet, 201

GetValvesActiveMap
 CWarnerValveControllerDeviceNet, 640

GetValvesManualStateMap
 CWarnerValveControllerDeviceNet, 640

GetValveTableEntry
 CWarnerValveControllerDeviceNet, 641

GetVDD3I
 CCMOSMea_FunctionNet, 105

GetVDDI
 CCMOSMea_FunctionNet, 105

GetVdsVgs
 CGrapheneFunctionNet, 169, 170

GetVdVs
 CGrapheneFunctionNet, 170

GetVdVsDAC
 CGrapheneFunctionNet, 170, 172

GetVersion
 CMcsUsbNet, 319

GetVersionInt
 DriverVersionNet, 672

GetVersionString
 DriverVersionNet, 673

GetVMMaxNegativeCurrent
 CMcsBus_VoltageModeNet, 233

GetVMMaxNegativeCurrentEeprom
 CMcsBus_VoltageModeNet, 233

GetVMMaxNegativeVoltage
 CMcsBus_VoltageModeNet, 233

GetVMMaxNegativeVoltageEeprom
 CMcsBus_VoltageModeNet, 233

GetVMMaxPositiveCurrent
 CMcsBus_VoltageModeNet, 233

GetVMMaxPositiveCurrentEeprom
 CMcsBus_VoltageModeNet, 233

GetVMMaxPositiveVoltage
 CMcsBus_VoltageModeNet, 233

GetVMMaxPositiveVoltageEeprom
 CMcsBus_VoltageModeNet, 234

GetVMOutputOnOff
 CMcsBus_VoltageModeNet, 234

GetVMVoltage
 CMcsBus_VoltageModeNet, 234

GetVolatileClampOffset
 CMultiwellDeviceNet, 378

GetVoltage
 CokuvisionStimulatorDeviceNet, 393

GetVoltage12V
 CRoboDeviceNet, 467

GetVoltage12VLimit
 CRoboDeviceNet, 467

GetVoltage5V
 CRoboDeviceNet, 467

GetVoltage5VLimit
 CRoboDeviceNet, 467

GetVoltageAirvalve
 CRoboDeviceNet, 467

GetVoltageAirvalveLimit
 CRoboDeviceNet, 467

GetVoltageClampControllerParam_D
 CWarnerUssingFunctionNet, 623

GetVoltageClampControllerParam_I
 CWarnerUssingFunctionNet, 624

GetVoltageClampControllerParam_P
 CWarnerUssingFunctionNet, 624

GetVoltageRange
 CGrapheneFunctionNet, 172

GetVoltageRangeIndex
 CMcsUsbDacqNet, 269

GetVoltageRangeInMicroVolt
 CMcsUsbDacqNet, 270

 CStg200xBasicNet, 525

 CStimulusFunctionNet, 561

 CW2100_StimulatorFunctionNet, 608

GetVoltageRangeInMilliVolt
 CMcsUsbDacqNet, 270

GetVoltageReached
 CGrapheneFunctionNet, 172, 173

GetVoltageResolution
 CGrapheneFunctionNet, 173

GetVoltageResolutionInMicroVolt
 CStg200xBasicNet, 525

 CStimulusFunctionNet, 562

 CW2100_StimulatorFunctionNet, 609

GetVoltageRs485A
 CRoboDeviceNet, 468

GetVoltageRs485ALimit
 CRoboDeviceNet, 468

GetVoltageRs485B
 CRoboDeviceNet, 468

GetVoltageRs485BLimit
 CRoboDeviceNet, 468

GetVoltageValves
 CRoboDeviceNet, 468

GetVoltageValvesLimit
 CRoboDeviceNet, 468

GetVolti
 CTcxDeviceNet, 580

GetWaveform
 CTEERFunctionNet, 589

GetWaveLengthInNanometer
 CMultiwellOptoStimFunctionNet, 384

GetWorkingFrequency
 CRFFFunctionNet, 442

GetWPADebugMode
 CWClassicFunctionNet, 660

GetWPAType
 CWClassicFunctionNet, 660

GetXGain

- CRoboDacqNet, [455](#)
- GetXilinxFlashOffset
 - CMcsUsbFactoryNet, [294](#)
- GetXilinxFlashReadCommand
 - CMcsUsbFactoryNet, [294](#)
- GND_SWITCH_BIT
 - CW2100_StimulatorFunctionNet, [611](#)
- Graphene_ASIC
 - Mcs::Usb, [75](#)
- GrapheneASIC
 - Mcs::Usb, [78](#)
- GrapheneASICHeadstage
 - Mcs::Usb, [65](#)
- GrapheneProjectTestDevice
 - Mcs::Usb, [76](#)
- Ground
 - Mcs::Usb, [54](#), [62](#)
- GroupID
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [121](#)
- GroupType
 - CDeviceGroupChannelInfoTemplateNet< Dacq-GroupChannelEnumTemplateNet >, [121](#)
- GyroOnly
 - Mcs::Usb, [87](#)
- Hardware
 - Mcs::Usb, [64](#)
- HasAnalogOut
 - CSCUFunctionNet, [503](#)
- HasGalvanicIsolation
 - CSCUFunctionNet, [503](#)
- HasHSPowerSwitch
 - CSCUFunctionNet, [503](#)
- HasIMU
 - HeadstageIDType, [678](#)
- HasOptoCurrentMeasurement
 - HeadstageIDType, [678](#)
- HasRadioControl
 - CRadioControlledDevicesNet, [435](#)
- HasRef
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- HasRefI
 - CRoboStatorDeviceNet, [485](#)
- HasRefXY
 - CRoboStatorDeviceNet, [485](#)
- HasRefZ
 - CRoboStatorDeviceNet, [485](#)
- HasSoftwareKey
 - CMcsUsbNet, [319](#)
- Headstage1NCBathCurrentGroup
 - Mcs::Usb, [57](#)
- Headstage1NCCol2CurrentGroup
 - Mcs::Usb, [57](#)
- Headstage1NChipTempGroup
 - Mcs::Usb, [57](#)
- HeadstageElectrodeGroup
 - Mcs::Usb, [57](#)
- HeadstageIDEnumNet
 - Mcs::Usb, [65](#)
- HeadstageIDType, [676](#)
 - CompareTo, [678](#)
 - ElectricalStimulation, [677](#)
 - Entry, [678](#)
 - Equals, [678](#)
 - HasIMU, [678](#)
 - HasOptoCurrentMeasurement, [678](#)
 - HeadstageIDType, [677](#)
 - HeadstageType, [678](#)
 - HeadstageTypeEnum, [677](#)
 - ID, [678](#)
 - MeasuringOnly, [677](#)
 - NumberOfAnalogChannels, [678](#)
 - NumberOfStimulationChannels, [679](#)
 - OpticalStimulation, [677](#)
 - SN, [679](#)
 - StimulusParameters, [679](#)
 - ToString, [678](#)
 - Type, [679](#)
 - TypeValue, [679](#)
 - Unknown, [677](#)
 - UserDefinedName, [679](#)
 - Valid, [679](#)
 - W16IsW14, [679](#)
- HeadstageIDTypeObject, [679](#)
 - _AdditionalText, [680](#)
 - _IdType, [680](#)
 - AdditionalText, [681](#)
 - Equals, [680](#)
 - GetHashCode, [680](#)
 - HeadstageIDTypeObject, [680](#)
 - IdType, [681](#)
 - ToString, [680](#)
- HeadstageIDTypeState, [681](#)
 - ControlState, [681](#)
 - DataState, [681](#)
 - IdType, [681](#)
 - State, [681](#)
- HeadstageType
 - HeadstageIDType, [678](#)
- HeadstageTypeEnum
 - HeadstageIDType, [677](#)
- HEKA_LIH3_DEVICE
 - Mcs::Usb, [59](#)
- HekaEPC10Double
 - Mcs::Usb, [74](#)
- HekaEPC10Quadro
 - Mcs::Usb, [74](#)
- HekaEPC10Single
 - Mcs::Usb, [74](#)
- HekaEPC10Triple
 - Mcs::Usb, [74](#)
- HekaEPCLite
 - Mcs::Usb, [74](#)
- HekaITEV100
 - Mcs::Usb, [74](#)

HekaLIH30
Mcs::Usb, [73](#)

HekaLIH406
Mcs::Usb, [74](#)

HekaLIH816
Mcs::Usb, [74](#)

HekaPG610
Mcs::Usb, [74](#)

HekaPG611
Mcs::Usb, [74](#)

HekaPG612
Mcs::Usb, [74](#)

HekaPG618
Mcs::Usb, [74](#)

HekaPG690
Mcs::Usb, [74](#)

HiClamp
Mcs::Usb, [76](#)

HiClamp4Uart
Mcs::Usb, [76](#)

Highpass
Mcs::Usb, [63](#)

HighSpeed
Mcs::Usb, [69](#)

HLA
Mcs::Usb, [74](#)

HLADacq
CHLADeviceNet, [179](#)

Hs1Digital
Mcs::Usb, [60](#)

HS1DigitalData1
Mcs::Usb, [83](#)

HS1ElectrodeGroup
Mcs::Usb, [69](#)

HS1Sideband1
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband10
Mcs::Usb, [70](#), [80](#)

HS1Sideband11
Mcs::Usb, [70](#), [80](#)

HS1Sideband12
Mcs::Usb, [70](#), [80](#)

HS1Sideband13
Mcs::Usb, [70](#)

HS1Sideband14
Mcs::Usb, [70](#)

HS1Sideband15
Mcs::Usb, [70](#)

HS1Sideband16
Mcs::Usb, [70](#)

HS1Sideband17
Mcs::Usb, [70](#)

HS1Sideband18
Mcs::Usb, [70](#)

HS1Sideband2
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband3
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband4
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband5
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband6
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Sideband7
Mcs::Usb, [70](#), [80](#)

HS1Sideband8
Mcs::Usb, [70](#), [80](#)

HS1Sideband9
Mcs::Usb, [70](#), [80](#)

Hs1SidebandHigh
Mcs::Usb, [60](#)

Hs1SidebandLow
Mcs::Usb, [60](#)

Hs1Trigger
Mcs::Usb, [60](#)

HS1Trigger10Status
Mcs::Usb, [70](#), [80](#)

HS1Trigger11Status
Mcs::Usb, [70](#), [80](#)

HS1Trigger12Status
Mcs::Usb, [70](#), [80](#)

HS1Trigger13Status
Mcs::Usb, [70](#)

HS1Trigger14Status
Mcs::Usb, [70](#)

HS1Trigger15Status
Mcs::Usb, [70](#)

HS1Trigger16Status
Mcs::Usb, [70](#)

HS1Trigger17Status
Mcs::Usb, [70](#)

HS1Trigger18Status
Mcs::Usb, [70](#)

HS1Trigger1Status
Mcs::Usb, [60](#), [70](#), [79](#)

HS1Trigger2Status
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Trigger3Status
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Trigger4Status
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Trigger5Status
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Trigger6Status
Mcs::Usb, [60](#), [70](#), [80](#)

HS1Trigger7Status
Mcs::Usb, [70](#), [80](#)

HS1Trigger8Status
Mcs::Usb, [70](#), [80](#)

HS1Trigger9Status
Mcs::Usb, [70](#), [80](#)

Hs2Digital
Mcs::Usb, [60](#)

HS2DigitalData1
Mcs::Usb, [83](#)

HS2ElectrodeGroup
 Mcs::Usb, [69](#)
 HS2Sideband1
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband10
 Mcs::Usb, [71](#), [80](#)
 HS2Sideband11
 Mcs::Usb, [71](#), [80](#)
 HS2Sideband12
 Mcs::Usb, [71](#), [80](#)
 HS2Sideband13
 Mcs::Usb, [71](#)
 HS2Sideband14
 Mcs::Usb, [71](#)
 HS2Sideband15
 Mcs::Usb, [71](#)
 HS2Sideband16
 Mcs::Usb, [71](#)
 HS2Sideband17
 Mcs::Usb, [71](#)
 HS2Sideband18
 Mcs::Usb, [71](#)
 HS2Sideband2
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband3
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband4
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband5
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband6
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Sideband7
 Mcs::Usb, [71](#), [80](#)
 HS2Sideband8
 Mcs::Usb, [71](#), [80](#)
 HS2Sideband9
 Mcs::Usb, [71](#), [80](#)
 Hs2SidebandHigh
 Mcs::Usb, [60](#)
 Hs2SidebandLow
 Mcs::Usb, [60](#)
 Hs2Trigger
 Mcs::Usb, [60](#)
 HS2Trigger10Status
 Mcs::Usb, [71](#), [80](#)
 HS2Trigger11Status
 Mcs::Usb, [71](#), [80](#)
 HS2Trigger12Status
 Mcs::Usb, [71](#), [80](#)
 HS2Trigger13Status
 Mcs::Usb, [71](#)
 HS2Trigger14Status
 Mcs::Usb, [71](#)
 HS2Trigger15Status
 Mcs::Usb, [71](#)
 HS2Trigger16Status
 Mcs::Usb, [71](#)
 HS2Trigger17Status
 Mcs::Usb, [71](#)
 HS2Trigger18Status
 Mcs::Usb, [71](#)
 HS2Trigger1Status
 Mcs::Usb, [60](#), [71](#), [80](#)
 HS2Trigger2Status
 Mcs::Usb, [60](#), [71](#), [80](#)
 HS2Trigger3Status
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Trigger4Status
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Trigger5Status
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Trigger6Status
 Mcs::Usb, [61](#), [71](#), [80](#)
 HS2Trigger7Status
 Mcs::Usb, [71](#), [80](#)
 HS2Trigger8Status
 Mcs::Usb, [71](#), [80](#)
 HS2Trigger9Status
 Mcs::Usb, [71](#), [80](#)
 HWInfo
 CMcsUsbDacqNet, [270](#)
 HwVersion
 CMcsUsbListEntryNet, [304](#)
 ID
 HeadStageIDType, [678](#)
 Idle
 Mcs::Usb, [81](#)
 IdProduct
 DeviceIdNet, [664](#)
 IdType
 HeadstageIDTypeObject, [681](#)
 HeadStageIDTypeState, [681](#)
 IdVendor
 DeviceIdNet, [664](#)
 IFB2GoldenInterfaceboard
 Mcs::Usb, [65](#)
 IFB30GoldenInterfaceboard
 Mcs::Usb, [65](#)
 IFChannel1
 Mcs::Usb, [56](#)
 IFChannel2
 Mcs::Usb, [56](#)
 IFChannel3
 Mcs::Usb, [56](#)
 IFChannel4
 Mcs::Usb, [56](#)
 IFChannel5
 Mcs::Usb, [56](#)
 IFChannel6
 Mcs::Usb, [56](#)
 IFChannel7
 Mcs::Usb, [56](#)
 IFChannel8
 Mcs::Usb, [56](#)
 IFDigChannelsGroup

- Mcs::Usb, [57](#), [69](#), [79](#), [89](#)
- IM16KRC
 - Mcs::Usb, [76](#)
- IM16S16KRA
 - Mcs::Usb, [76](#)
- IM16S8KRA
 - Mcs::Usb, [76](#)
- IM64KRB
 - Mcs::Usb, [76](#)
- IM64KRC
 - Mcs::Usb, [76](#)
- Input
 - Mcs::Usb, [73](#)
- Intel
 - Mcs::Usb, [86](#)
- InterfaceADCGroup
 - Mcs::Usb, [57](#), [69](#), [79](#), [88](#)
- InterfaceBoard2
 - Mcs::Usb, [65](#)
- IntToDouble
 - Mcs::Usb, [63](#)
- InvitroSignalCollectorUnit
 - Mcs::Usb, [65](#)
- InvivoSignalCollectorUnit
 - Mcs::Usb, [65](#)
- IoVoltageEnumNet
 - Mcs::Usb, [65](#)
- IS32KRA
 - Mcs::Usb, [76](#)
- IsAnalogOutEnabled
 - CSCUFunctionNet, [503](#)
- IsAutomaticAnalogOut
 - CSCUFunctionNet, [503](#)
- IsBusy
 - CPPCFunctionNet, [419](#)
- IsChamberAvailable
 - CWarnerUssingFunctionNet, [624](#)
- IsChipPowered
 - CCMOSMea_FunctionNet, [105](#)
- IsConnected
 - CMcsUsbNet, [319](#)
- IsDeviceHighSpeed
 - CMcsUsbNet, [319](#)
- IsDeviceHighSpeedCapable
 - CMcsUsbNet, [319](#)
- IsDeviceTypeOf
 - CMcsUsbListNet, [307](#)
- IsDigitalChannelDedicated
 - CMcsUsbDacqNet::CHWInfo, [182](#)
- IsDigitalOutPortInverted
 - CWarnerValveControllerDeviceNet, [641](#)
- IsDigitalOutPortInvertedEvent
 - CWarnerValveControllerDeviceNet, [653](#)
- IsEqual
 - CFilterCoefficientsNet, [130](#)
- IsExceptionsEnabled
 - CMcsUsbNet, [320](#)
- IsGateFloating
 - CCMOSMea_FunctionNet, [105](#)
- IsHeadstageAvailable
 - CSCUFunctionNet, [503](#)
- IsHeadstageAvailableEvent
 - CSCUFunctionNet, [507](#)
- IsHighCurrentMode
 - CWarnerUssingFunctionNet, [625](#)
- IsHSPowered
 - CSCUFunctionNet, [504](#)
- IsInDacqLegacyMode
 - CSCUFunctionNet, [504](#)
- IsInternalCalibrationFinished
 - CTEERFunctionNet, [590](#)
 - CWarnerUssingFunctionNet, [625](#)
- IsPlateTypeValid
 - CMultiwellDeviceNet, [379](#)
- IsPulseEnabled
 - CWarnerUssingFunctionNet, [625](#)
- IsPumpMotorOn
 - CRoboFluidDeviceNet, [479](#)
- IsQueueEnabled
 - CRoboDeviceNet, [468](#)
- IsQueueStarted
 - CRoboDeviceNet, [468](#)
- IsRunning
 - CMeaCleanDeviceNet, [336](#)
 - CMeaCoatDeviceNet, [341](#)
- IsSamplingFinished
 - CTEERFunctionNet, [590](#)
- IsUserTriggerEnabled
 - CLIH3DeviceNet, [193](#)
- IsValveDigitalInInverted
 - CWarnerValveControllerDeviceNet, [641](#)
- IsValveDigitalInInvertedEvent
 - CWarnerValveControllerDeviceNet, [653](#)
- IsValveOpen
 - CWarnerValveControllerDeviceNet, [642](#)
- IsValveOpenEvent
 - CWarnerValveControllerDeviceNet, [654](#)
- IsValveOpenInAnalogMode
 - CWarnerValveControllerDeviceNet, [642](#)
- IsValveOpenInAnalogModeEvent
 - CWarnerValveControllerDeviceNet, [654](#)
- IsValveOpenInDigitalMode
 - CWarnerValveControllerDeviceNet, [642](#)
- IsValveOpenInDigitalModeEvent
 - CWarnerValveControllerDeviceNet, [654](#)
- Kelvin
 - Mcs::Usb, [51](#)
- LastPosition
 - Mcs::Usb, [61](#), [71](#), [80](#), [83](#), [89](#)
- LIH30_ADC_Channel_EnumNet
 - Mcs::Usb, [66](#)
- LIH30_DAC_Channel_EnumNet
 - Mcs::Usb, [66](#)
- LIH30_EPC10_Bus_EnumNet
 - Mcs::Usb, [66](#)

- LIH30ADCCtrl
 - Mcs::Usb, [65](#)
- LIH30ADCModulesGroup
 - Mcs::Usb, [57](#)
- LIH30Interfaceboard
 - Mcs::Usb, [65](#)
- LIH30TestADCGroup
 - Mcs::Usb, [57](#)
- LIH30UserADCGroup
 - Mcs::Usb, [57](#)
- ListModeSendStart
 - CStg200xBasicNet, [526](#)
- ListModeSendStop
 - CStg200xBasicNet, [526](#)
- ListOfChangedTriggers
 - StgStatusNet, [695](#)
- LoadPressure
 - CPPCFunctionNet, [421](#)
- LoadUserFirmware
 - CMcsUsbFactoryNet, [294](#), [295](#)
- LoadValveTable
 - CWarnerValveControllerDeviceNet, [642](#)
- Lock
 - Mcs::Usb, [73](#)
- LockPlateClamp
 - CMultiwellDeviceNet, [379](#)
- Lowpass
 - Mcs::Usb, [63](#)
- LowSpeed
 - Mcs::Usb, [69](#)
- m_Bottom
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [436](#)
- m_Left
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [436](#)
- m_pGilsonDevice
 - CGilsonDeviceNet, [164](#)
- m_pMcsBus_MotorControlNet
 - CRoboFluidDeviceNet, [480](#)
- m_pMcsUsb
 - CMcsUsbFunctionNet, [299](#)
- m_pMcsUsbFunction
 - CMcsUsbFunctionNet, [299](#)
- m_pRoboFluidDevice
 - CRoboFluidDeviceNet, [480](#)
- m_Right
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [436](#)
- m_Top
 - CCMOSMeaDeviceNet::CRegionOfInterestRect, [436](#)
- Manual
 - Mcs::Usb, [73](#), [81](#), [90](#)
- Manufacturer
 - CMcsUsbListEntryNet, [304](#)
- MaxBitNumber
 - DigitalSource< digitalsourceenum >, [665](#)
- DigitalSourceGeneral, [666](#), [667](#)
- MaxBitNumberStatic
 - DigitalSource< digitalsourceenum >, [665](#)
- MBC08
 - Mcs::Usb, [75](#)
- MbcChannelStateEnumNet
 - Mcs::Usb, [67](#)
- MbcChargingModeEnumNet
 - Mcs::Usb, [67](#)
- MbcRatedCapacityEnumNet
 - Mcs::Usb, [67](#)
- MC_Card
 - Mcs::Usb, [73](#)
- MCS
 - Mcs::Usb, [87](#)
- Mcs, [22](#)
- Mcs::Usb, [22](#)
 - A, [67](#)
 - AccelOnly, [87](#)
 - AdapterTypeEnumNet, [50](#)
 - ALA_VC3, [73](#), [87](#)
 - ALA_VC3_DEVICE, [59](#)
 - ALTERA, [52](#)
 - ALTERA_BASE, [53](#)
 - ALTERA_BOOTSTRAP, [53](#)
 - ALTERA_GOLD, [53](#)
 - ALTERA_TARGET1, [54](#)
 - ALTERA_TARGET2, [54](#)
 - ALTERA_TARGET3, [54](#)
 - AlwaysOn, [76](#)
 - Ampere, [51](#), [87](#)
 - AMS_Dongle, [74](#)
 - Analog, [73](#), [90](#)
 - AnalogGroup, [58](#)
 - AnalogOut_DAC_Range_EnumNet, [51](#)
 - AnalogSource_HS1, [51](#)
 - AnalogSource_HS2, [51](#)
 - AnalogSource_IF, [51](#)
 - AnalogSourceEnumNet, [51](#)
 - AnalogUnitEnumNet, [51](#)
 - Any, [73](#), [87](#)
 - Armed, [81](#)
 - ASMedia, [86](#)
 - AudioTestChannelGroup, [57](#), [69](#), [79](#), [89](#)
 - Aux, [60](#)
 - AuxIn, [60](#), [70](#), [79](#), [82](#), [89](#)
 - AuxPort, [56](#)
 - B, [67](#)
 - Bessel, [64](#)
 - BMI, [83](#)
 - Bootstrap, [52](#)
 - BootstrapOtherCypress, [52](#)
 - Both, [87](#)
 - Break, [77](#)
 - BUS1MCSBUS1, [52](#)
 - BUS1MCSBUS2, [52](#)
 - BUS2MCSBUS1, [52](#)
 - BUS2MCSBUS2, [52](#)

- BUSNUMBER1, [52](#)
- BUSNUMBER2, [52](#)
- Butterworth, [64](#)
- Campden_Ci4600EphysVideoDataIntegrator, [73](#)
- CatchAmp, [72](#)
- CFirmwareDestinationNet, [51](#)
- channeldata_current, [82](#)
- channeldata_current_always_boost, [82](#)
- channeldata_current_always_boost_own_sync, [82](#)
- channeldata_current_own_boost_gnd_sync, [82](#)
- channeldata_current_own_sync, [82](#)
- channeldata_positive_current, [82](#)
- channeldata_positive_current_own_boost_gnd_sync, [82](#)
- channeldata_positive_current_own_sync, [82](#)
- channeldata_positive_voltage, [82](#)
- channeldata_voltage, [82](#)
- ChannelPIC, [52](#)
- ChannelTest, [75](#)
- ChecksumAndPacketCounter, [56](#)
- Ci4600Intan, [50](#)
- ClampModeCurrent, [84](#)
- ClampModeInternalCalibration, [84](#)
- ClampModeOpen, [84](#)
- ClampModeVoltage, [83](#)
- Close, [72](#)
- CmosMea, [65](#)
- CMOSMeaBathModeEnumNet, [54](#)
- CmosMeaHeadstage, [64](#)
- CMOSMeaHeadstage1NCBathCurrentEnumNet, [54](#)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet, [54](#)
- CMOSMeaHeadstage1NChipTempEnumNet, [55](#)
- CMOSMeaHS1SidebandEnumNet, [55](#)
- CMOSMeaHS1TriggerStatusEnumNet, [55](#)
- CmosmealFB2, [65](#)
- CMOSMealFDigChannelEnumNet, [55](#)
- CMOSMealInterfaceADCEnumNet, [56](#)
- CmosMealInterfaceboard, [64](#)
- CMOSMeaPacketFrameContextGroupEnumNet, [56](#)
- CMOSMeaSTG1DACSignalEnumNet, [56](#)
- CMOSMeaValueUnitEnumNet, [57](#)
- CommaPositionA, [63](#)
- CommaPositionB, [63](#)
- csCapacityTestDischarge, [67](#)
- csCapacityTestPrecharge, [67](#)
- csCharge, [67](#)
- csDischarge, [67](#)
- csError, [67](#)
- csIdleChargeFinished, [67](#)
- csIdleNoBattery, [67](#)
- csRefreshBattery, [67](#)
- CurrentClamp, [87](#)
- CurrentMeasure, [54](#)
- Cypress, [87](#)
- Cypress_FX1, [73](#)
- Cypress_FX2, [73](#)
- Cypress_FX3, [73](#)
- DAC1Channel, [56](#)
- DAC2Channel, [56](#)
- DAC3Channel, [56](#)
- DAC4Channel, [56](#)
- DACQ1DigitalGroup, [57](#)
- DacqGroupChannelEnumNet, [57](#)
- DacqMeaGroupTypeEnumNet, [57](#)
- DacqTrigger, [62](#)
- DataModeEnumNet, [58](#)
- DEST_FX3_TARGET_MASK, [54](#)
- DEST_TARGET1, [53](#)
- DEST_TARGET10, [53](#)
- DEST_TARGET11, [53](#)
- DEST_TARGET12, [53](#)
- DEST_TARGET13, [54](#)
- DEST_TARGET14, [54](#)
- DEST_TARGET15, [54](#)
- DEST_TARGET2, [53](#)
- DEST_TARGET3, [53](#)
- DEST_TARGET4, [53](#)
- DEST_TARGET5, [53](#)
- DEST_TARGET6, [53](#)
- DEST_TARGET7, [53](#)
- DEST_TARGET8, [53](#)
- DEST_TARGET9, [53](#)
- DEST_TARGET_MASK, [54](#)
- DEVICE_NOT_FOUND, [63](#)
- DeviceEnumNet, [58](#)
- DeviceHasNoHeadstage, [65, 78](#)
- DeviceNotConnected, [64, 65, 78](#)
- DeviceRunStatus, [61, 70, 79, 83, 89](#)
- DigDataFromReceiver, [89](#)
- Digital, [73, 90](#)
- DigitalData, [61, 70, 79, 83, 89](#)
- DigitalDataStreamEnableEnumNet, [59](#)
- DigitalGroup, [58](#)
- DigitalIn, [59, 60, 70, 79, 82, 89](#)
- DigitalInOfOutPort, [60, 70, 79, 82, 89](#)
- DigitalInPort, [55](#)
- DigitalInReserverd, [59](#)
- DigitalMux, [55](#)
- DigitalOut, [59](#)
- DigitalOutReg, [55](#)
- DigitalOutReserved, [60](#)
- DigitalOutStimulator, [61, 70, 79, 83, 89](#)
- DigitalPulse, [60, 70, 79, 82, 89](#)
- DigitalReg, [56](#)
- DigitalSourceEnumNet, [60](#)
- DigitalStimulatorTriggerEventEnumNet, [61](#)
- DigitalStimulatorTriggerSlopeEnumNet, [61](#)
- DigitalTargetEnumNet, [61](#)
- Digout, [62](#)
- DigOutStim, [60](#)
- DigOutStimulatorStartTrigger, [62](#)
- DigOutStimulatorStopTrigger, [62](#)
- Digstream, [62](#)

- DigStreamFromReceiver, [89](#)
- DigStreamToReceiver, [62](#)
- Dilutor, [76](#)
- DongleS, [74](#)
- Dotriapot, [74](#)
- DoubleToInt, [63](#)
- DownloadOnly, [81](#)
- DSP, [51](#)
- DSPAnalogGroup, [69](#), [79](#)
- DSPDataGroup, [57](#), [88](#)
- DSPDigitalGroup, [69](#), [79](#)
- eCube, [75](#), [78](#)
- eCubeHeadstage, [65](#)
- ElectrodeDacMuxEnumNet, [62](#)
- ElectrodeModeEnumNet, [62](#)
- ElectrodeOffset, [87](#)
- emAutomatic, [62](#)
- emManual, [62](#)
- Encapsulator, [76](#)
- enCMosMeaChipType, [62](#)
- EnSTG200x_STATUS, [63](#)
- EOFAndCRC, [56](#)
- ExternBCTester, [74](#)
- ExternDTester, [74](#)
- ExternSTester, [74](#)
- Falling, [61](#)
- FCB, [74](#)
- FCX, [74](#)
- Feedback, [60](#), [70](#), [79](#), [82](#), [89](#)
- FeedbackHigh, [60](#)
- FeedbackLow, [60](#)
- FeedbackReg, [56](#)
- FilterAttributeEnumNet, [63](#)
- FilterBandEnumNet, [63](#)
- FilterCalculationDirectionEnumNet, [63](#)
- FilterFamilyEnumNet, [63](#)
- FilterTypeEnumNet, [64](#)
- Finished, [81](#)
- FPGA10, [52](#)
- FPGA10_BASE, [53](#)
- FPGA10_GOLD, [53](#)
- FPGA11, [52](#)
- FPGA11_BASE, [53](#)
- FPGA11_GOLD, [53](#)
- FPGA12, [52](#)
- FPGA12_BASE, [53](#)
- FPGA12_GOLD, [53](#)
- FPGA13, [52](#)
- FPGA13_BASE, [53](#)
- FPGA13_GOLD, [53](#)
- FPGA14, [52](#)
- FPGA14_BASE, [53](#)
- FPGA14_GOLD, [53](#)
- FPGA15, [52](#)
- FPGA15_BASE, [53](#)
- FPGA15_GOLD, [53](#)
- FPGA16, [53](#)
- FPGA16_BASE, [53](#)
- FPGA16_GOLD, [53](#)
- FPGA2, [52](#)
- FPGA2_BASE, [53](#)
- FPGA2_GOLD, [53](#)
- FPGA3, [52](#)
- FPGA3_BASE, [53](#)
- FPGA3_GOLD, [53](#)
- FPGA4, [52](#)
- FPGA4_BASE, [53](#)
- FPGA4_GOLD, [53](#)
- FPGA5, [52](#)
- FPGA5_BASE, [53](#)
- FPGA5_GOLD, [53](#)
- FPGA6, [52](#)
- FPGA6_BASE, [53](#)
- FPGA6_GOLD, [53](#)
- FPGA7, [52](#)
- FPGA7_BASE, [53](#)
- FPGA7_GOLD, [53](#)
- FPGA8, [52](#)
- FPGA8_BASE, [53](#)
- FPGA8_GOLD, [53](#)
- FPGA9, [52](#)
- FPGA9_BASE, [53](#)
- FPGA9_GOLD, [53](#)
- FPGA_BASE, [53](#)
- FPGA_BOOTSTRAP, [53](#)
- FPGA_GOLD, [53](#)
- FPGA_NORMAL, [51](#)
- FpgaldEnumNet, [64](#)
- FrameContextGroup, [58](#)
- FullCharge, [67](#)
- FullSpeed, [69](#)
- FunkDongleS, [74](#)
- Gated_High_Active, [76](#)
- Gated_Low_Active, [76](#)
- GE2100, [75](#)
- Generic, [74](#)
- Graphene_ASIC, [75](#)
- GrapheneASIC, [78](#)
- GrapheneASICHeadstage, [65](#)
- GrapheneProjectTestDevice, [76](#)
- Ground, [54](#), [62](#)
- GyroOnly, [87](#)
- Hardware, [64](#)
- Headstage1NCBathCurrentGroup, [57](#)
- Headstage1NCCol2CurrentGroup, [57](#)
- Headstage1NChipTempGroup, [57](#)
- HeadstageElectrodeGroup, [57](#)
- HeadstageIdEnumNet, [65](#)
- HEKA_LIH3_DEVICE, [59](#)
- HekaEPC10Double, [74](#)
- HekaEPC10Quadro, [74](#)
- HekaEPC10Single, [74](#)
- HekaEPC10Triple, [74](#)
- HekaEPCLite, [74](#)
- HekaITEV100, [74](#)
- HekaLIH30, [73](#)

HekaLIH406, [74](#)
HekaLIH816, [74](#)
HekaPG610, [74](#)
HekaPG611, [74](#)
HekaPG612, [74](#)
HekaPG618, [74](#)
HekaPG690, [74](#)
HiClamp, [76](#)
HiClamp4Uart, [76](#)
Highpass, [63](#)
HighSpeed, [69](#)
HLA, [74](#)
Hs1Digital, [60](#)
Hs1DigitalData1, [83](#)
Hs1ElectrodeGroup, [69](#)
Hs1Sideband1, [60](#), [70](#), [80](#)
Hs1Sideband10, [70](#), [80](#)
Hs1Sideband11, [70](#), [80](#)
Hs1Sideband12, [70](#), [80](#)
Hs1Sideband13, [70](#)
Hs1Sideband14, [70](#)
Hs1Sideband15, [70](#)
Hs1Sideband16, [70](#)
Hs1Sideband17, [70](#)
Hs1Sideband18, [70](#)
Hs1Sideband2, [60](#), [70](#), [80](#)
Hs1Sideband3, [60](#), [70](#), [80](#)
Hs1Sideband4, [60](#), [70](#), [80](#)
Hs1Sideband5, [60](#), [70](#), [80](#)
Hs1Sideband6, [60](#), [70](#), [80](#)
Hs1Sideband7, [70](#), [80](#)
Hs1Sideband8, [70](#), [80](#)
Hs1Sideband9, [70](#), [80](#)
Hs1SidebandHigh, [60](#)
Hs1SidebandLow, [60](#)
Hs1Trigger, [60](#)
Hs1Trigger10Status, [70](#), [80](#)
Hs1Trigger11Status, [70](#), [80](#)
Hs1Trigger12Status, [70](#), [80](#)
Hs1Trigger13Status, [70](#)
Hs1Trigger14Status, [70](#)
Hs1Trigger15Status, [70](#)
Hs1Trigger16Status, [70](#)
Hs1Trigger17Status, [70](#)
Hs1Trigger18Status, [70](#)
Hs1Trigger1Status, [60](#), [70](#), [79](#)
Hs1Trigger2Status, [60](#), [70](#), [80](#)
Hs1Trigger3Status, [60](#), [70](#), [80](#)
Hs1Trigger4Status, [60](#), [70](#), [80](#)
Hs1Trigger5Status, [60](#), [70](#), [80](#)
Hs1Trigger6Status, [60](#), [70](#), [80](#)
Hs1Trigger7Status, [70](#), [80](#)
Hs1Trigger8Status, [70](#), [80](#)
Hs1Trigger9Status, [70](#), [80](#)
Hs2Digital, [60](#)
Hs2DigitalData1, [83](#)
Hs2ElectrodeGroup, [69](#)
Hs2Sideband1, [61](#), [71](#), [80](#)
Hs2Sideband10, [71](#), [80](#)
Hs2Sideband11, [71](#), [80](#)
Hs2Sideband12, [71](#), [80](#)
Hs2Sideband13, [71](#)
Hs2Sideband14, [71](#)
Hs2Sideband15, [71](#)
Hs2Sideband16, [71](#)
Hs2Sideband17, [71](#)
Hs2Sideband18, [71](#)
Hs2Sideband2, [61](#), [71](#), [80](#)
Hs2Sideband3, [61](#), [71](#), [80](#)
Hs2Sideband4, [61](#), [71](#), [80](#)
Hs2Sideband5, [61](#), [71](#), [80](#)
Hs2Sideband6, [61](#), [71](#), [80](#)
Hs2Sideband7, [71](#), [80](#)
Hs2Sideband8, [71](#), [80](#)
Hs2Sideband9, [71](#), [80](#)
Hs2SidebandHigh, [60](#)
Hs2SidebandLow, [60](#)
Hs2Trigger, [60](#)
Hs2Trigger10Status, [71](#), [80](#)
Hs2Trigger11Status, [71](#), [80](#)
Hs2Trigger12Status, [71](#), [80](#)
Hs2Trigger13Status, [71](#)
Hs2Trigger14Status, [71](#)
Hs2Trigger15Status, [71](#)
Hs2Trigger16Status, [71](#)
Hs2Trigger17Status, [71](#)
Hs2Trigger18Status, [71](#)
Hs2Trigger1Status, [60](#), [71](#), [80](#)
Hs2Trigger2Status, [60](#), [71](#), [80](#)
Hs2Trigger3Status, [61](#), [71](#), [80](#)
Hs2Trigger4Status, [61](#), [71](#), [80](#)
Hs2Trigger5Status, [61](#), [71](#), [80](#)
Hs2Trigger6Status, [61](#), [71](#), [80](#)
Hs2Trigger7Status, [71](#), [80](#)
Hs2Trigger8Status, [71](#), [80](#)
Hs2Trigger9Status, [71](#), [80](#)
Idle, [81](#)
IFB2GoldenInterfaceboard, [65](#)
IFB30GoldenInterfaceboard, [65](#)
IFChannel1, [56](#)
IFChannel2, [56](#)
IFChannel3, [56](#)
IFChannel4, [56](#)
IFChannel5, [56](#)
IFChannel6, [56](#)
IFChannel7, [56](#)
IFChannel8, [56](#)
IFDigChannelsGroup, [57](#), [69](#), [79](#), [89](#)
IM16KRC, [76](#)
IM16S16KRA, [76](#)
IM16S8KRA, [76](#)
IM64KRB, [76](#)
IM64KRC, [76](#)
Input, [73](#)
Intel, [86](#)
InterfaceADCGroup, [57](#), [69](#), [79](#), [88](#)

InterfaceBoard2, 65
 IntToDouble, 63
 InvitroSignalCollectorUnit, 65
 InvivoSignalCollectorUnit, 65
 IoVoltageEnumNet, 65
 IS32KRA, 76
 Kelvin, 51
 LastPosition, 61, 71, 80, 83, 89
 LIH30_ADC_Channel_EnumNet, 66
 LIH30_DAC_Channel_EnumNet, 66
 LIH30_EPC10_Bus_EnumNet, 66
 LIH30ADCCtrl, 65
 LIH30ADCModulesGroup, 57
 LIH30Interfaceboard, 65
 LIH30TestADCGroup, 57
 LIH30UserADCGroup, 57
 Lock, 73
 Lowpass, 63
 LowSpeed, 69
 Manual, 73, 81, 90
 MBC08, 75
 MbcChannelStateEnumNet, 67
 MbcChargingModeEnumNet, 67
 MbcRatedCapacityEnumNet, 67
 MC_Card, 73
 MCS, 87
 MCS_ANY_BUS, 69
 MCS_CHANNELTEST_DEVICE, 58
 MCS_DEVICE_ANY, 58
 MCS_DEVICE_USB, 58
 MCS_DEVICE_USB_CYPRESS, 59
 MCS_ENCAPSULATOR_DEVICE, 59
 MCS_EXTERN_BC_TESTER_DEVICE, 59
 MCS_EXTERN_D_TESTER_DEVICE, 59
 MCS_FCX_DEVICE, 58
 MCS_FYI_DEVICE, 59
 MCS_GENERIC_DEVELOPMENT_DEVICE, 58
 MCS_HICLAMP_DEVICE, 59
 MCS_HLA_DEVICE, 59
 MCS_MBC08_DEVICE, 59
 MCS_MC_STIMULUS_DEVICE, 58
 MCS_MCCARD_DEVICE, 58
 MCS_MEA_CLEAN_DEVICE, 59
 MCS_MEA_COAT_DEVICE, 59
 MCS_MEA_DEVICE, 58
 MCS_MEA_IMPEDANCE_DEVICE, 58
 MCS_MEA_SWITCH_DEVICE, 58
 MCS_MEASURETABLE_DEVICE, 59
 MCS_MEAUSB_DEVICE, 58
 MCS_NF_GEN_DEVICE, 59
 MCS_OCTOPOT_DEVICE, 58
 MCS_OKUVISION_STIMULATOR_DEVICE, 59
 MCS_PATCHSERVER_DEVICE, 59
 MCS_PATHIDENT_DEVICE, 59
 MCS_PCI_BUS, 69
 MCS_PCX_DEVICE, 58
 MCS_PEDOTER_DEVICE, 59
 MCS_PERISTALTIC_PUMP_DEVICE, 59
 MCS_PGA_DEVICE, 58
 MCS_PPC_DEVICE, 59
 MCS_PPS5_DEVICE, 59
 MCS_PPS_DEVICE, 59
 MCS_RETINA_AMS_DONGLE, 59
 MCS_RETINA_LED_DEVICE, 58
 MCS_ROBO_DEVICE, 59
 MCS_ROBOINJECT_DEVICE, 59
 MCS_ROBOOCYTE2_DEVICE, 59
 MCS_SAFEIS_DEVICE, 59
 MCS_SMARTIMPLANT_DEVICE, 59
 MCS_SOFTWARE_DONGLE_DEVICE, 59
 MCS_STG_DEVICE, 58
 MCS_SW2TO64_DEVICE, 58
 MCS_TCX_DEVICE, 58
 MCS_TERSENS_DEVICE, 58
 MCS_UNDEFINED_BUS, 69
 MCS_USB_BUS, 69
 MCSBUS0, 52
 MCSBUS1, 51
 MCSBUS10, 52
 MCSBUS11, 52
 MCSBUS12, 52
 MCSBUS13, 52
 MCSBUS14, 52
 MCSBUS15, 52
 MCSBUS2, 51
 MCSBUS3, 52
 MCSBUS4, 52
 MCSBUS5, 52
 MCSBUS6, 52
 MCSBUS7, 52
 MCSBUS8, 52
 MCSBUS9, 52
 McsBusTypeEnumNet, 67
 McsUsbSpeedEnumNet, 69
 MCU1, 51
 ME128, 75
 ME16, 75
 ME2100, 75
 Me2100_32PICiCE40, 78
 Me2100_32PICiCE40Headstage, 64
 Me2100_32Xilinx, 78
 Me2100_32XilinxHeadstage, 64
 Me2100Graphene16_32, 78
 Me2100Graphene16_32Headstage, 65
 Me2100Interfaceboard, 64
 Me2100InvitroSignalCollectorUnit, 64
 Me2100InvivoSignalCollectorUnit, 64
 Me2100UPA32, 78
 Me2100UPA32Headstage, 64
 ME256, 75
 ME32, 75
 ME64, 75
 MEA1060, 75
 MEA120, 50
 MEA2100, 75
 Mea2100, 65

MEA2100_256, 75
Mea2100_256, 65
MEA2100_256DacqGroupChannelEnumNet, 69
MEA2100_256DigitalSourceEnumNet, 69
Mea2100_256Headstage, 64
Mea2100_256Interfaceboard, 64
MEA2100_32, 75
MEA2100_Lite, 75
Mea2100_Lite, 65
MEA2100_Mini, 75
MEA2100_Mini_Usb_develop, 75
MEA2100BetaScreen, 75
Mea2100BetaScreen, 78
Mea2100BetaScreenHeadstage, 64
Mea2100Headstage, 64
Mea2100Interfaceboard, 64
Mea2100LiteHeadstage, 65
Mea2100Mini120, 78
Mea2100Mini120Headstage, 64
Mea2100Mini60ECP5, 78
Mea2100Mini60ECP5Headstage, 65
Mea2100Mini60PICIcE40, 78
Mea2100Mini60PICIcE40Headstage, 64
Mea2100MultiwellIFB2, 64
Mea2100STG, 64
MEA252, 50
MEA256, 75
MEA2x32, 50
MEA2x60, 50
MEA32, 50
MEA60, 50
MEA_2_252_2, 50
MEA_2_252_2_6Well, 50
MEA_2_252_2_9Well, 50
MEA_2_252_2_Test, 50
MEA_Clean, 75
MEA_Coat, 75
MEA_Impedance, 75
MEA_Sanofi, 75
MEA_Switch, 75
MEA_Switch_2_1, 75
MEA_Switch_4_2, 75
MeaLayoutEnumNet, 71
MeasureTable, 76
MicroAmpere, 57
MilliDegreeCelsius, 57
mIMEA60, 71
mIUnknown, 71
ModulA_ADC0, 66
ModulA_ADC1, 66
ModulA_ADC2, 66
ModulA_ADC3, 66
ModulA_DAC0, 66
ModulA_DAC1, 66
ModulB_ADC0, 66
ModulB_ADC1, 66
ModulB_ADC2, 66
ModulB_ADC3, 66
ModulB_DAC0, 66
ModulB_DAC1, 66
ModulC_ADC0, 66
ModulC_ADC1, 66
ModulC_ADC2, 66
ModulC_ADC3, 66
ModulC_DAC0, 66
ModulC_DAC1, 66
ModulD_ADC0, 66
ModulD_ADC1, 66
ModulD_ADC2, 66
ModulD_ADC3, 66
ModulD_DAC0, 66
ModulD_DAC1, 66
Monitor, 81
Movement, 77
Multiboot, 75
Multiwell, 65, 75
Multiwell96, 50
Multiwell_ICC, 75
Multiwell_MEA_Mini, 75
MultiwellHeadstage, 64
MultiwellInterfaceboard, 64
MultiwellMini, 78
MultiwellMiniHeadstage, 64
MultiwellOptoStim, 74
MultiwellPlateTypeEnumNet, 72
Mux, 59
MuxOtherDevice, 59
Nanion, 83
NanoAmpere, 57
NanoVolt, 57
NCBathCurrent, 54
NCCol2Current, 54
NChipTemperature, 55
Neptun, 76
NeuroChip, 75
NeurochipConfig, 76
NF_Gen, 75
NineWell, 82
nMos16LV, 62
nMos32LV, 62
nMos36LN, 62
nMos64LN, 62
No_Plate, 72
None, 50, 59, 73, 81, 86, 87
Normal, 72
NOT_CONNECTED, 63
NotApplicable, 50
NoUnit, 57
NTC10K, 83
Octopot, 74
Off, 76, 77, 87
off, 77
OK, 63
Okuvision_Stimulator, 74
OnChannelData, 90
OnDeviceArrivalRemoval, 90

- One, [60](#), [70](#), [79](#), [82](#), [89](#)
- OnError, [90](#)
- OnMcsUsbDeviceState, [90](#)
- OnMcsUsbDeviceStateCallback, [91](#)
- OnMwPollStatus, [91](#)
- OnStg200xDataHandler, [91](#)
- OnStg200xErrorHandler, [91](#)
- OnStgPollStatus, [91](#)
- OnUpdateFirmwareProgress, [91](#)
- OnUpdateFirmwareStatusChange, [91](#)
- Open, [72](#)
- OpenClamp, [87](#)
- Output, [73](#)
- PacketFrameContextGroup, [57](#), [69](#), [79](#), [89](#)
- PatchServAdcModeEnumNet, [72](#)
- PatchServer, [76](#)
- PathIdent, [75](#)
- PC, [90](#)
- PCI, [87](#)
- PCX, [74](#)
- PeriodicPulse, [60](#)
- PeristalticPump, [75](#)
- PGA, [74](#)
- PIC, [52](#)
- PIC10, [52](#)
- PIC11, [52](#)
- PIC12, [52](#)
- PIC2, [52](#)
- PIC3, [52](#)
- PIC4, [52](#)
- PIC5, [52](#)
- PIC6, [52](#)
- PIC7, [52](#)
- PIC8, [52](#)
- PIC9, [52](#)
- PicoAmpere, [57](#)
- Plate_24W030MGA, [72](#)
- Plate_24W300_30_1152GBA, [72](#)
- Plate_24W300_30GBA, [72](#)
- Plate_24W300_30GBB, [72](#)
- Plate_24W300_30GMA, [72](#)
- Plate_24W700_100FMA, [72](#)
- Plate_24W700_100FMB, [72](#)
- Plate_24W700_100FMC, [72](#)
- Plate_24W700_100PBA, [72](#)
- Plate_72W500_100FMA, [72](#)
- Plate_72W500_100PMA, [72](#)
- Plate_96W300_80_1152FMA, [72](#)
- Plate_96W400_80_1152FMB, [72](#)
- Plate_96W700_100FMA, [72](#)
- Plate_96W700_100FMB, [72](#)
- Plate_96W700_100GBA, [72](#)
- Plate_96W700_100GBB, [72](#)
- Plate_96W700_100GBC, [72](#)
- Plate_96W700_100GBD, [72](#)
- Plate_96W700_100GMA, [72](#)
- Plate_Dummy, [72](#)
- Plate_Dummy_126, [72](#)
- PlateClampEnumNet, [72](#)
- PlateClampLockEnumNet, [73](#)
- PlusMinus10Volts, [51](#)
- PlusMinus2Comma5Volts, [51](#)
- PlusMinus5Volts, [51](#)
- PortDirectionEnumNet, [73](#)
- Pos900, [76](#)
- PositionBase, [76](#)
- PositionIIBase, [76](#)
- PositionIICentralUnit, [76](#)
- PositionImp, [76](#)
- PostCommaA, [63](#)
- PostCommaB, [63](#)
- PP_Pump_Mode_Type_EnumNet, [73](#)
- PPC, [75](#)
- PPS2, [75](#)
- PPS4plus1, [75](#)
- PPS5, [75](#)
- PPS5_DIG, [75](#)
- PreCommaA, [63](#)
- PreCommaB, [63](#)
- ProductIdEnumNet, [73](#)
- PT100, [83](#)
- PT1000, [83](#)
- PulseGenerator, [61](#), [70](#), [79](#), [83](#), [89](#)
- PulseGenerator_Mode_EnumNet, [76](#)
- raGate, [77](#)
- raIgnore, [77](#)
- raRestart, [77](#)
- raSingle, [77](#)
- raStop, [77](#)
- rawdata, [82](#)
- RC, [64](#)
- rc100mAh, [67](#)
- rc200mAh, [67](#)
- rc300mAh, [67](#)
- rc30mAh, [67](#)
- rcGreater300mAh, [67](#)
- Rectangle, [84](#)
- Ref16, [77](#)
- Ref24, [77](#)
- Ref32, [77](#)
- Ref8, [77](#)
- Reference, [77](#)
- ReferenceElectrodeModeEnumNet, [76](#)
- ReferenceElectrodeSwitchPositionEnumNet, [77](#)
- RegisterHigh, [60](#)
- RegisterLow, [60](#)
- Regular, [83](#)
- Renesas, [86](#)
- Reserved1, [83](#)
- Reserved2, [83](#)
- Reserved3, [83](#)
- Reserved4, [83](#)
- Reserved5, [83](#)
- Retina_LED, [74](#)
- RetriggerActionEnumNet, [77](#)
- Rising, [61](#)

RoboCurrentModeEnumNet, [77](#)
RoboInject, [76](#)
Roboocyte2, [76](#)
RoboStatusEventDelegate, [91](#)
Running, [81](#)
SafeIS, [76](#)
SampleDstSize16, [78](#)
SampleDstSize32, [78](#)
SampleDstSizeNet, [78](#)
SampleSize16Signed, [78](#)
SampleSize16Unsigned, [78](#)
SampleSize24Signed, [78](#)
SampleSize24Unsigned, [78](#)
SampleSize32Signed, [78](#)
SampleSize32Unsigned, [78](#)
SampleSize64Signed, [78](#)
SampleSize64Unsigned, [78](#)
SampleSizeNet, [78](#)
SBSVector1, [55](#)
SBSVector2, [55](#)
SBSVector3, [55](#)
SBSVector4, [55](#)
SCU1ElectrodeGroupHS1, [79](#)
SCU1ElectrodeGroupHS2, [79](#)
SCU1ElectrodeGroupHS3, [79](#)
SCU1ElectrodeGroupHS4, [79](#)
SCU2ElectrodeGroupHS1, [79](#)
SCU2ElectrodeGroupHS2, [79](#)
SCU2ElectrodeGroupHS3, [79](#)
SCU2ElectrodeGroupHS4, [79](#)
SCU_HeadstageIdEnumNet, [78](#)
SCUDacqGroupChannelEnumNet, [79](#)
SCUDigitalSourceEnumNet, [79](#)
Settings, [90](#)
Signed_16bit, [58](#)
Signed_24bit, [58](#)
Signed_32bit, [58](#)
Sine, [84](#)
SingleWell, [82](#)
SixWell, [82](#)
SmartImplant, [76](#)
SOFAAndCTRLword, [56](#)
Software, [64](#)
SoftwareDongle, [75](#)
Standby, [77](#), [87](#)
Start, [61](#)
State, [87](#)
STG, [74](#)
Stg1, [62](#)
STG1DACSignalGroup, [57](#), [69](#), [79](#)
STG1SidebandsGroup, [57](#), [69](#), [79](#)
STG1TriggerStatusGroup, [57](#), [69](#), [79](#)
Stg2, [62](#)
Stg200xDigoutModeEnumNet, [81](#)
Stg200xSegmentFlagsEnumNet, [81](#)
Stg200xTriggerStatusEnumNet, [81](#)
STG2DACSignalGroup, [69](#), [79](#)
STG2SidebandsGroup, [69](#), [79](#)
STG2TriggerStatusGroup, [69](#), [79](#)
Stg3, [62](#)
STG3008_FA, [74](#)
STG4002, [74](#)
STG4002_opto, [74](#)
STG4004, [74](#)
STG4004_opto, [74](#)
STG4008, [74](#)
STG4008_opto, [74](#)
STG400x, [74](#)
STG400x_opto, [74](#)
STG5, [74](#)
STG_DestinationEnumNet, [82](#)
StgListModeTrigger, [62](#)
StgTrigger, [62](#)
Stimulation, [54](#)
StimulationLayoutConfigurationEnumNet, [82](#)
Stop, [61](#), [72](#)
StorageCharge, [67](#)
SubtractFromAll, [77](#)
SubtractFromAllOther, [77](#)
SubtractFromReferenceElectrodeOnly, [77](#)
SubtractionOff, [77](#)
SuperSpeed, [69](#)
Sw2to64, [75](#)
SYNCOUT1, [81](#)
SYNCOUT2, [81](#)
SYNCOUT3, [81](#)
SYNCOUT4, [81](#)
SYNCOUT5, [81](#)
SYNCOUT6, [81](#)
SYNCOUT7, [81](#)
SYNCOUT8, [81](#)
syncoutdata, [82](#)
SyncStart, [81](#)
Table, [90](#)
TBSI_127, [50](#)
TBSI_15, [50](#)
TBSI_31, [50](#)
TBSI_5, [50](#)
TBSI_63, [50](#)
TBSI_Dacq, [75](#)
TBSI_DACQDigitalSourceEnumNet, [82](#)
TBSI_Reserved, [50](#)
TbsiDacq, [65](#)
TbsiDacqHeadstage, [64](#)
TbsiDacqInterfaceboard, [64](#)
TC01, [74](#)
TC02, [74](#)
TCX, [74](#)
TcxDeviceTypeEnumNet, [83](#)
TcxSensorTypeEnumNet, [83](#)
TeerClampModeEnumNet, [83](#)
TeerWaveformEnumNet, [84](#)
Tersens, [74](#)
Test_ADC_EPC10, [66](#)
Test_DAC_EPC10, [66](#)
Timestamp, [56](#)

TouchTest, [90](#)
Triggerbox_AMS, [74](#)
Triggerbox_AMS3, [74](#)
Triggerbox_IMS, [74](#)
Triggerbox_R5, [74](#)
TriggerOnly, [81](#)
TriggerSourceEnumNet, [84](#)
TriggerStatus1, [55](#)
TriggerStatus2, [55](#)
TriggerStatus3, [55](#)
TriggerStatus4, [55](#)
tsAuxIn1, [85](#)
tsAuxIn2, [85](#)
tsDACQCy1Dev1Runs, [86](#)
tsDACQCy1Dev2Runs, [86](#)
tsDACQCy2Dev1Runs, [86](#)
tsDACQCy2Dev2Runs, [86](#)
tsDigitalIn1, [84](#)
tsDigitalIn10, [84](#)
tsDigitalIn11, [84](#)
tsDigitalIn12, [84](#)
tsDigitalIn13, [84](#)
tsDigitalIn14, [84](#)
tsDigitalIn15, [84](#)
tsDigitalIn16, [84](#)
tsDigitalIn17, [84](#)
tsDigitalIn18, [84](#)
tsDigitalIn19, [84](#)
tsDigitalIn2, [84](#)
tsDigitalIn20, [84](#)
tsDigitalIn21, [84](#)
tsDigitalIn22, [84](#)
tsDigitalIn23, [84](#)
tsDigitalIn24, [84](#)
tsDigitalIn25, [85](#)
tsDigitalIn26, [85](#)
tsDigitalIn27, [85](#)
tsDigitalIn28, [85](#)
tsDigitalIn29, [85](#)
tsDigitalIn3, [84](#)
tsDigitalIn30, [85](#)
tsDigitalIn31, [85](#)
tsDigitalIn32, [85](#)
tsDigitalIn4, [84](#)
tsDigitalIn5, [84](#)
tsDigitalIn6, [84](#)
tsDigitalIn7, [84](#)
tsDigitalIn8, [84](#)
tsDigitalIn9, [84](#)
tsDigitalPuse0, [85](#)
tsDigitalPuse1, [85](#)
tsDigitalPuse10, [86](#)
tsDigitalPuse11, [86](#)
tsDigitalPuse12, [86](#)
tsDigitalPuse13, [86](#)
tsDigitalPuse14, [86](#)
tsDigitalPuse15, [86](#)
tsDigitalPuse16, [86](#)
tsDigitalPuse17, [86](#)
tsDigitalPuse18, [86](#)
tsDigitalPuse19, [86](#)
tsDigitalPuse2, [85](#)
tsDigitalPuse20, [86](#)
tsDigitalPuse21, [86](#)
tsDigitalPuse22, [86](#)
tsDigitalPuse23, [86](#)
tsDigitalPuse24, [86](#)
tsDigitalPuse25, [86](#)
tsDigitalPuse26, [86](#)
tsDigitalPuse27, [86](#)
tsDigitalPuse28, [86](#)
tsDigitalPuse29, [86](#)
tsDigitalPuse3, [85](#)
tsDigitalPuse30, [86](#)
tsDigitalPuse31, [86](#)
tsDigitalPuse4, [85](#)
tsDigitalPuse5, [85](#)
tsDigitalPuse6, [85](#)
tsDigitalPuse7, [85](#)
tsDigitalPuse8, [86](#)
tsDigitalPuse9, [86](#)
tsFeedback1, [85](#)
tsFeedback10, [85](#)
tsFeedback11, [85](#)
tsFeedback12, [85](#)
tsFeedback13, [85](#)
tsFeedback14, [85](#)
tsFeedback15, [85](#)
tsFeedback16, [85](#)
tsFeedback17, [85](#)
tsFeedback18, [85](#)
tsFeedback19, [85](#)
tsFeedback2, [85](#)
tsFeedback20, [85](#)
tsFeedback21, [85](#)
tsFeedback22, [85](#)
tsFeedback23, [85](#)
tsFeedback24, [85](#)
tsFeedback25, [85](#)
tsFeedback26, [85](#)
tsFeedback27, [85](#)
tsFeedback28, [85](#)
tsFeedback29, [85](#)
tsFeedback3, [85](#)
tsFeedback30, [85](#)
tsFeedback31, [85](#)
tsFeedback32, [85](#)
tsFeedback4, [85](#)
tsFeedback5, [85](#)
tsFeedback6, [85](#)
tsFeedback7, [85](#)
tsFeedback8, [85](#)
tsFeedback9, [85](#)
tsNone, [84](#)
tsSidebandBit8, [86](#)
tsTriggered, [86](#)

Unknown, 50, 51, 63, 64, 83, 86
unknown, 62
UnknownDest, 54
UnknownSpeed, 69
Unlock, 73
Unsigned_16bit, 58
Unsigned_24bit, 58
Unsigned_32bit, 58
UpdateTrigger, 81
USB, 51
USB_TARGET1, 54
USB_TARGET2, 54
USB_TARGET3, 54
UsbTest, 75
UsbVendorIdEnumNet, 86
User_ADC_0, 66
User_ADC_1, 66
User_ADC_2, 66
User_ADC_3, 66
User_ADC_4, 66
User_DAC_0, 66
User_DAC_1, 66
User_DAC_2, 66
UssingChamber, 65
UssingClampModeEnumNet, 86
UssingRail, 65
UssingUnitEnumNet, 87
VendorIdEnumNet, 87
Volt, 51, 87
Voltage_3V3, 65
Voltage_5V0, 65
VoltageClamp, 87
W2100, 75
W2100_Accel_Gyro_Select_EnumNet, 87
W2100DacqGroupChannelEnumNet, 87
W2100DigitalSourceEnumNet, 89
W2100IFB2, 65
W2100Interfaceboard, 64
W2100WirelessReceiver, 64, 65
W2100WirelessReceiverAnalog, 64, 65
Warner, 83
Warner_TEER_Machine, 76
Warner_Ussing, 76
WARNER_USSING_DEVICE, 59
Warner_Valve_Control, 76
WARNER_VALVE_CONTROL_DEVICE, 59
Whole_Cell_Patch, 75
WholeCellPatch, 78
WholeCellPatchHeadstage, 65
WirelessHeadStageAccDataRE1HS1, 88
WirelessHeadStageAccDataRE1HS2, 88
WirelessHeadStageAccDataRE1HS3, 88
WirelessHeadStageAccDataRE1HS4, 88
WirelessHeadStageAccDataRE2HS1, 88
WirelessHeadStageAccDataRE2HS2, 88
WirelessHeadStageAccDataRE2HS3, 88
WirelessHeadStageAccDataRE2HS4, 88
WirelessHeadStageAnalogRE1HS1, 88
WirelessHeadStageAnalogRE1HS2, 88
WirelessHeadStageAnalogRE1HS3, 88
WirelessHeadStageAnalogRE1HS4, 88
WirelessHeadStageAnalogRE2HS1, 88
WirelessHeadStageAnalogRE2HS2, 88
WirelessHeadStageAnalogRE2HS3, 88
WirelessHeadStageAnalogRE2HS4, 88
WirelessHeadStageGyroDataRE1HS1, 88
WirelessHeadStageGyroDataRE1HS2, 88
WirelessHeadStageGyroDataRE1HS3, 88
WirelessHeadStageGyroDataRE1HS4, 88
WirelessHeadStageGyroDataRE2HS1, 88
WirelessHeadStageGyroDataRE2HS2, 88
WirelessHeadStageGyroDataRE2HS3, 88
WirelessHeadStageGyroDataRE2HS4, 88
WirelessHeadStageOptoStimCurrentRE1HS1, 88
WirelessHeadStageOptoStimCurrentRE1HS2, 88
WirelessHeadStageOptoStimCurrentRE1HS3, 88
WirelessHeadStageOptoStimCurrentRE1HS4, 88
WirelessHeadStageOptoStimCurrentRE2HS1, 88
WirelessHeadStageOptoStimCurrentRE2HS2, 89
WirelessHeadStageOptoStimCurrentRE2HS3, 89
WirelessHeadStageOptoStimCurrentRE2HS4, 89
WirelessHeadStageReservedARE1HS1, 88
WirelessHeadStageReservedARE1HS2, 88
WirelessHeadStageReservedARE1HS3, 88
WirelessHeadStageReservedARE1HS4, 88
WirelessHeadStageReservedARE2HS1, 89
WirelessHeadStageReservedARE2HS2, 89
WirelessHeadStageReservedARE2HS3, 89
WirelessHeadStageReservedARE2HS4, 89
WirelessHeadStageReservedBRE1HS1, 89
WirelessHeadStageReservedBRE1HS2, 89
WirelessHeadStageReservedBRE1HS3, 89
WirelessHeadStageReservedBRE1HS4, 89
WirelessHeadStageReservedBRE2HS1, 89
WirelessHeadStageReservedBRE2HS2, 89
WirelessHeadStageReservedBRE2HS3, 89
WirelessHeadStageReservedBRE2HS4, 89
WirelessHeadStageReservedCRE1HS1, 89
WirelessHeadStageReservedCRE1HS2, 89
WirelessHeadStageReservedCRE1HS3, 89
WirelessHeadStageReservedCRE1HS4, 89
WirelessHeadStageReservedCRE2HS1, 89
WirelessHeadStageReservedCRE2HS2, 89
WirelessHeadStageReservedCRE2HS3, 89
WirelessHeadStageReservedCRE2HS4, 89
WirelessHeadStageStatusRE1HS1, 88
WirelessHeadStageStatusRE1HS2, 88
WirelessHeadStageStatusRE1HS3, 88
WirelessHeadStageStatusRE1HS4, 88
WirelessHeadStageStatusRE2HS1, 88
WirelessHeadStageStatusRE2HS2, 88
WirelessHeadStageStatusRE2HS3, 88
WirelessHeadStageStatusRE2HS4, 88
WirelessTestAdapter, 50
Work, 90
WPA16, 75

- WPA32, [75](#)
- WPA4, [75](#)
- WPA8, [75](#)
- WvcDisplayModeEnumNet, [90](#)
- WvcValveModeEnumNet, [90](#)
- Zero, [60](#), [70](#), [79](#), [82](#), [89](#)
- MCS_ANY_BUS
 - Mcs::Usb, [69](#)
- MCS_CHANNELTEST_DEVICE
 - Mcs::Usb, [58](#)
- MCS_DEVICE_ANY
 - Mcs::Usb, [58](#)
- MCS_DEVICE_USB
 - Mcs::Usb, [58](#)
- MCS_DEVICE_USB_CYPRESS
 - Mcs::Usb, [59](#)
- MCS_ENCAPSULATOR_DEVICE
 - Mcs::Usb, [59](#)
- MCS_EXTERN_BC_TESTER_DEVICE
 - Mcs::Usb, [59](#)
- MCS_EXTERN_D_TESTER_DEVICE
 - Mcs::Usb, [59](#)
- MCS_FCX_DEVICE
 - Mcs::Usb, [58](#)
- MCS_FYI_DEVICE
 - Mcs::Usb, [59](#)
- MCS_GENERIC_DEVELOPMENT_DEVICE
 - Mcs::Usb, [58](#)
- MCS_HICLAMP_DEVICE
 - Mcs::Usb, [59](#)
- MCS_HLA_DEVICE
 - Mcs::Usb, [59](#)
- MCS_MBC08_DEVICE
 - Mcs::Usb, [59](#)
- MCS_MC_STIMULUS_DEVICE
 - Mcs::Usb, [58](#)
- MCS_MCCARD_DEVICE
 - Mcs::Usb, [58](#)
- MCS_MEA_CLEAN_DEVICE
 - Mcs::Usb, [59](#)
- MCS_MEA_COAT_DEVICE
 - Mcs::Usb, [59](#)
- MCS_MEA_DEVICE
 - Mcs::Usb, [58](#)
- MCS_MEA_IMPEDANCE_DEVICE
 - Mcs::Usb, [58](#)
- MCS_MEA_SWITCH_DEVICE
 - Mcs::Usb, [58](#)
- MCS_MEASURETABLE_DEVICE
 - Mcs::Usb, [59](#)
- MCS_MEAUSB_DEVICE
 - Mcs::Usb, [58](#)
- MCS_NF_GEN_DEVICE
 - Mcs::Usb, [59](#)
- MCS_OCTOPOT_DEVICE
 - Mcs::Usb, [58](#)
- MCS_OKUVISION_STIMULATOR_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PATCHSERVER_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PATHIDENT_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PCI_BUS
 - Mcs::Usb, [69](#)
- MCS_PCX_DEVICE
 - Mcs::Usb, [58](#)
- MCS_PEDOTER_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PERISTALTIC_PUMP_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PGA_DEVICE
 - Mcs::Usb, [58](#)
- MCS_PPC_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PPS5_DEVICE
 - Mcs::Usb, [59](#)
- MCS_PPS_DEVICE
 - Mcs::Usb, [59](#)
- MCS_RETINA_AMS_DONGLE
 - Mcs::Usb, [59](#)
- MCS_RETINA_LED_DEVICE
 - Mcs::Usb, [58](#)
- MCS_ROBO_DEVICE
 - Mcs::Usb, [59](#)
- MCS_ROBOINJECT_DEVICE
 - Mcs::Usb, [59](#)
- MCS_ROBOOCYTE2_DEVICE
 - Mcs::Usb, [59](#)
- MCS_SAFEIS_DEVICE
 - Mcs::Usb, [59](#)
- MCS_SMARTIMPLANT_DEVICE
 - Mcs::Usb, [59](#)
- MCS_SOFTWARE_DONGLE_DEVICE
 - Mcs::Usb, [59](#)
- MCS_STG_DEVICE
 - Mcs::Usb, [58](#)
- MCS_SW2TO64_DEVICE
 - Mcs::Usb, [58](#)
- MCS_TCX_DEVICE
 - Mcs::Usb, [58](#)
- MCS_TERSENS_DEVICE
 - Mcs::Usb, [58](#)
- MCS_UNDEFINED_BUS
 - Mcs::Usb, [69](#)
- MCS_USB_BUS
 - Mcs::Usb, [69](#)
- McsBus
 - CPPCDeviceNet, [414](#)
 - CPPS_DeviceNet, [423](#)
 - CRoboDeviceNet, [476](#)
- MCSBUS0
 - Mcs::Usb, [52](#)
- MCSBUS1
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [51](#)
- MCSBUS10

- FirmwareDestinationNames, [675](#)
- Mcs::Usb, [52](#)
- MCSBUS11
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS12
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS13
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS14
 - Mcs::Usb, [52](#)
- MCSBUS15
 - Mcs::Usb, [52](#)
- MCSBUS2
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [51](#)
- MCSBUS3
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS4
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS5
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS6
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS7
 - FirmwareDestinationNames, [675](#)
 - Mcs::Usb, [52](#)
- MCSBUS8
 - FirmwareDestinationNames, [676](#)
 - Mcs::Usb, [52](#)
- MCSBUS9
 - FirmwareDestinationNames, [676](#)
 - Mcs::Usb, [52](#)
- McsBus_MotorControl
 - CPeristalticPumpDeviceNet, [399](#)
 - CPPCDeviceNet, [414](#)
 - CPPS_DeviceNet, [423](#)
 - CRoboDeviceNet, [476](#)
 - CRoboFluidDeviceNet, [480](#)
- McsBus_Sensor
 - CPPCDeviceNet, [414](#)
 - CPPS_DeviceNet, [423](#)
- McsBus_VoltageMode
 - CFluidControlDeviceNet, [143](#)
- McsBus_XY
 - CRoboDeviceNet, [473](#)
- McsBus_ZI
 - CRoboDeviceNet, [473](#)
- McsBusTypeEnumNet
 - Mcs::Usb, [67](#)
- McsUsbDeviceStateEvent
 - CMcsUsbDeviceStatePushFunctionNet, [288](#)
 - CMcsUsbDeviceStatePushNet, [289](#)
- McsUsbSpeedEnumNet
 - Mcs::Usb, [69](#)
- MCU1
 - FirmwareDestinationNames, [676](#)
 - Mcs::Usb, [51](#)
- ME128
 - Mcs::Usb, [75](#)
- ME16
 - Mcs::Usb, [75](#)
- ME2100
 - Mcs::Usb, [75](#)
- Me2100_32PICI40
 - Mcs::Usb, [78](#)
- Me2100_32PICI40Headstage
 - Mcs::Usb, [64](#)
- Me2100_32Xilinx
 - Mcs::Usb, [78](#)
- Me2100_32XilinxHeadstage
 - Mcs::Usb, [64](#)
- Me2100Graphene16_32
 - Mcs::Usb, [78](#)
- Me2100Graphene16_32Headstage
 - Mcs::Usb, [65](#)
- Me2100Interfaceboard
 - Mcs::Usb, [64](#)
- Me2100InvitroSignalCollectorUnit
 - Mcs::Usb, [64](#)
- Me2100InvivoSignalCollectorUnit
 - Mcs::Usb, [64](#)
- Me2100UPA32
 - Mcs::Usb, [78](#)
- Me2100UPA32Headstage
 - Mcs::Usb, [64](#)
- ME256
 - Mcs::Usb, [75](#)
- ME32
 - Mcs::Usb, [75](#)
- ME64
 - Mcs::Usb, [75](#)
- MEA1060
 - Mcs::Usb, [75](#)
- MEA120
 - Mcs::Usb, [50](#)
- MEA2100
 - Mcs::Usb, [75](#)
- Mea2100
 - Mcs::Usb, [65](#)
- MEA2100_256
 - Mcs::Usb, [75](#)
- Mea2100_256
 - Mcs::Usb, [65](#)
- MEA2100_256DacqGroupChannelEnumNet
 - Mcs::Usb, [69](#)
- MEA2100_256DigitalSourceEnumNet
 - Mcs::Usb, [69](#)
- Mea2100_256Headstage
 - Mcs::Usb, [64](#)

- Mea2100_256Interfaceboard
 - Mcs::Usb, [64](#)
- MEA2100_32
 - Mcs::Usb, [75](#)
- MEA2100_Lite
 - Mcs::Usb, [75](#)
- Mea2100_Lite
 - Mcs::Usb, [65](#)
- MEA2100_Mini
 - Mcs::Usb, [75](#)
- MEA2100_Mini_Usb_develop
 - Mcs::Usb, [75](#)
- MEA2100BetaScreen
 - Mcs::Usb, [75](#)
- Mea2100BetaScreen
 - Mcs::Usb, [78](#)
- Mea2100BetaScreenHeadstage
 - Mcs::Usb, [64](#)
- Mea2100Headstage
 - Mcs::Usb, [64](#)
- Mea2100Interfaceboard
 - Mcs::Usb, [64](#)
- Mea2100LiteHeadstage
 - Mcs::Usb, [65](#)
- Mea2100Mini120
 - Mcs::Usb, [78](#)
- Mea2100Mini120Headstage
 - Mcs::Usb, [64](#)
- Mea2100Mini60ECP5
 - Mcs::Usb, [78](#)
- Mea2100Mini60ECP5Headstage
 - Mcs::Usb, [65](#)
- Mea2100Mini60PICiGE40
 - Mcs::Usb, [78](#)
- Mea2100Mini60PICiGE40Headstage
 - Mcs::Usb, [64](#)
- Mea2100MultiwellIFB2
 - Mcs::Usb, [64](#)
- Mea2100STG
 - Mcs::Usb, [64](#)
- MEA252
 - Mcs::Usb, [50](#)
- MEA256
 - Mcs::Usb, [75](#)
- MEA2x32
 - Mcs::Usb, [50](#)
- MEA2x60
 - Mcs::Usb, [50](#)
- MEA32
 - Mcs::Usb, [50](#)
- MEA60
 - Mcs::Usb, [50](#)
- MEA_2_252_2
 - Mcs::Usb, [50](#)
- MEA_2_252_2_6Well
 - Mcs::Usb, [50](#)
- MEA_2_252_2_9Well
 - Mcs::Usb, [50](#)
- MEA_2_252_2_Test
 - Mcs::Usb, [50](#)
- MEA_Clean
 - Mcs::Usb, [75](#)
- MEA_Coat
 - Mcs::Usb, [75](#)
- MEA_Impedance
 - Mcs::Usb, [75](#)
- MEA_Sanofi
 - Mcs::Usb, [75](#)
- MEA_Switch
 - Mcs::Usb, [75](#)
- MEA_Switch_2_1
 - Mcs::Usb, [75](#)
- MEA_Switch_4_2
 - Mcs::Usb, [75](#)
- MeaAudioFunctionNet
 - CMeaDeviceNet, [351](#)
- MeaDigitalDataFunctionNet
 - CMeaDeviceNet, [351](#)
- MeaFeedbackFunctionNet
 - CMeaDeviceNet, [351](#)
- MeaLayoutEnumNet
 - Mcs::Usb, [71](#)
- Measure
 - CPathIdentDeviceNet, [396](#)
- MeasureReservoir
 - CPPCFunctionNet, [421](#)
- MeasureTable
 - Mcs::Usb, [76](#)
- MeasuringOnly
 - HeadStageIDType, [677](#)
- MeFunctionNet
 - CMeaDeviceNet, [351](#)
- MicroAmpere
 - Mcs::Usb, [57](#)
- MilliDegreeCelsius
 - Mcs::Usb, [57](#)
- mkfilter
 - mkfilterNet, [682](#)
- mkfilter_coef_in_one_set
 - mkfilterNet, [682](#)
- mkfilter_highpass_coeff
 - mkfilterNet, [683](#)
- mkfilter_highpass_frequency_from_coeff
 - mkfilterNet, [683](#)
- mkfilter_highpass_frequency_from_k
 - mkfilterNet, [683](#)
- mkfilter_highpass_k
 - mkfilterNet, [683](#)
- mkfilter_MCS
 - mkfilterNet, [683](#)
- mkfilter_MCS_k
 - mkfilterNet, [683](#), [684](#)
- mkfilter_normalize_coeffs_int
 - mkfilterNet, [684](#)
- mkfilter_normalize_coeffs_short
 - mkfilterNet, [684](#)

- mkfilter_normalize_scale_coeffs_int
 - mkfilterNet, [684](#)
- mkfilter_scale_coef_in_one_set
 - mkfilterNet, [684](#)
- mkfilterNet, [682](#)
 - mkfilter, [682](#)
 - mkfilter_coef_in_one_set, [682](#)
 - mkfilter_highpass_coeff, [683](#)
 - mkfilter_highpass_frequency_from_coeff, [683](#)
 - mkfilter_highpass_frequency_from_k, [683](#)
 - mkfilter_highpass_k, [683](#)
 - mkfilter_MCS, [683](#)
 - mkfilter_MCS_k, [683](#), [684](#)
 - mkfilter_normalize_coeffs_int, [684](#)
 - mkfilter_normalize_coeffs_short, [684](#)
 - mkfilter_normalize_scale_coeffs_int, [684](#)
 - mkfilter_scale_coef_in_one_set, [684](#)
- mlMEA60
 - Mcs::Usb, [71](#)
- mlUnknown
 - Mcs::Usb, [71](#)
- ModulA_ADC0
 - Mcs::Usb, [66](#)
- ModulA_ADC1
 - Mcs::Usb, [66](#)
- ModulA_ADC2
 - Mcs::Usb, [66](#)
- ModulA_ADC3
 - Mcs::Usb, [66](#)
- ModulA_DAC0
 - Mcs::Usb, [66](#)
- ModulA_DAC1
 - Mcs::Usb, [66](#)
- ModulB_ADC0
 - Mcs::Usb, [66](#)
- ModulB_ADC1
 - Mcs::Usb, [66](#)
- ModulB_ADC2
 - Mcs::Usb, [66](#)
- ModulB_ADC3
 - Mcs::Usb, [66](#)
- ModulB_DAC0
 - Mcs::Usb, [66](#)
- ModulB_DAC1
 - Mcs::Usb, [66](#)
- ModulC_ADC0
 - Mcs::Usb, [66](#)
- ModulC_ADC1
 - Mcs::Usb, [66](#)
- ModulC_ADC2
 - Mcs::Usb, [66](#)
- ModulC_ADC3
 - Mcs::Usb, [66](#)
- ModulC_DAC0
 - Mcs::Usb, [66](#)
- ModulC_DAC1
 - Mcs::Usb, [66](#)
- ModulD_ADC0
 - Mcs::Usb, [66](#)
- ModulD_ADC1
 - Mcs::Usb, [66](#)
- ModulD_ADC2
 - Mcs::Usb, [66](#)
- ModulD_ADC3
 - Mcs::Usb, [66](#)
- ModulD_DAC0
 - Mcs::Usb, [66](#)
- ModulD_DAC1
 - Mcs::Usb, [66](#)
- Monitor
 - Mcs::Usb, [81](#)
- MoveAbs
 - CRoboDeviceNet, [468](#), [469](#)
- MoveAbsI
 - CRoboStatorDeviceNet, [485](#)
- MoveAbsXY
 - CRoboStatorDeviceNet, [485](#)
- MoveAbsZ
 - CRoboStatorDeviceNet, [485](#), [486](#)
- Movement
 - Mcs::Usb, [77](#)
- Multiboot
 - Mcs::Usb, [75](#)
- MultibootGetCypressImageId
 - CMcsUsbNet, [320](#)
- MultibootGetImageId
 - CMcsUsbNet, [320](#)
- MultibootGetSelectedImage
 - CMcsUsbNet, [320](#)
- MultibootSelectImage
 - CMcsUsbNet, [320](#)
- Multiwell
 - Mcs::Usb, [65](#), [75](#)
- MultiwellI96
 - Mcs::Usb, [50](#)
- Multiwell_ICC
 - Mcs::Usb, [75](#)
- Multiwell_MEA_Mini
 - Mcs::Usb, [75](#)
- MultiwellHeadstage
 - Mcs::Usb, [64](#)
- MultiwellInterfaceboard
 - Mcs::Usb, [64](#)
- MultiwellMini
 - Mcs::Usb, [78](#)
- MultiwellMiniHeadstage
 - Mcs::Usb, [64](#)
- MultiwellOptoStim
 - Mcs::Usb, [74](#)
- MultiwellPlateTypeEnumNet
 - Mcs::Usb, [72](#)
- Mux
 - Mcs::Usb, [59](#)
- MuxOtherDevice
 - Mcs::Usb, [59](#)
- MwPollStatusEvent

- CStg200xDownloadNet, 555
- Nanion
 - Mcs::Usb, 83
- NanoAmpere
 - Mcs::Usb, 57
- NanoVolt
 - Mcs::Usb, 57
- NCBathCurrent
 - Mcs::Usb, 54
- NCCol2Current
 - Mcs::Usb, 54
- NChipTemperature
 - Mcs::Usb, 55
- Neptun
 - Mcs::Usb, 76
- NeuroChip
 - Mcs::Usb, 75
- NeurochipConfig
 - Mcs::Usb, 76
- NF_Gen
 - Mcs::Usb, 75
- NineWell
 - Mcs::Usb, 82
- nMos16LV
 - Mcs::Usb, 62
- nMos32LV
 - Mcs::Usb, 62
- nMos36LN
 - Mcs::Usb, 62
- nMos64LN
 - Mcs::Usb, 62
- No_Plate
 - Mcs::Usb, 72
- None
 - Mcs::Usb, 50, 59, 73, 81, 86, 87
- Normal
 - Mcs::Usb, 72
- NOT_CONNECTED
 - Mcs::Usb, 63
- NotApplicable
 - Mcs::Usb, 50
- NoUnit
 - Mcs::Usb, 57
- NTC10K
 - Mcs::Usb, 83
- NullCommand
 - CRoboDeviceNet, 470
- NumberOfAnalogChannels
 - HeadStageIDType, 678
- NumberOfChannels
 - CDeviceGroupChannelInfoTemplateNet < Dacq-GroupChannelEnumTemplateNet >, 121
- NumberOfStimulationChannels
 - HeadStageIDType, 679
- NumCoefSets
 - CCreateFilterNet, 112
- Octopot
 - Mcs::Usb, 74
- Off
 - Mcs::Usb, 76, 77, 87
- off
 - Mcs::Usb, 77
- OK
 - Mcs::Usb, 63
- Okuvision_Stimulator
 - Mcs::Usb, 74
- OnChannelData
 - Mcs::Usb, 90
- OnDeviceArrivalRemoval
 - Mcs::Usb, 90
- One
 - Mcs::Usb, 60, 70, 79, 82, 89
- OnError
 - Mcs::Usb, 90
- OnGetActiveRunningTableNumber
 - CWarnerValveControllerDeviceNet, 643
- OnGetAnalogThresholdHigh
 - CWarnerValveControllerDeviceNet, 643
- OnGetAnalogThresholdLow
 - CWarnerValveControllerDeviceNet, 643
- OnGetAnalogVoltage
 - CWarnerValveControllerDeviceNet, 643
- OnGetAvailableHeadstages
 - CSCUFunctionNet, 504
- OnGetCurrentNumberOfValves
 - CWarnerValveControllerDeviceNet, 643
- OnGetDigitalOutPortValve
 - CWarnerValveControllerDeviceNet, 643
- OnGetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, 643
- OnGetDisplayMode
 - CWarnerValveControllerDeviceNet, 643
- OnGetPlateClampStateByHeadstage
 - CMultiwellCallbackFunctionNet, 374
- OnGetTableNamebyIndex
 - CWarnerValveControllerDeviceNet, 643
- OnGetValveActive
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveBoardRevision
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveLedOn
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveManualGroup
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveManualState
 - CWarnerValveControllerDeviceNet, 644
- OnGetValveMode
 - CWarnerValveControllerDeviceNet, 644
- OnIsDigitalOutPortInverted
 - CWarnerValveControllerDeviceNet, 644
- OnIsHeadstageAvailable
 - CSCUFunctionNet, 504
- OnIsValveDigitalInInverted

- CWarnerValveControllerDeviceNet, 645
- OnIsValveOpen
 - CWarnerValveControllerDeviceNet, 645
- OnIsValveOpenInAnalogMode
 - CWarnerValveControllerDeviceNet, 645
- OnIsValveOpenInDigitalMode
 - CWarnerValveControllerDeviceNet, 645
- OnMcsUsbDeviceState
 - Mcs::Usb, 90
- OnMcsUsbDeviceStateCallback
 - Mcs::Usb, 91
- OnMwPollStatus
 - Mcs::Usb, 91
- OnStg200xDataHandler
 - Mcs::Usb, 91
- OnStg200xErrorHandler
 - Mcs::Usb, 91
- OnStgPollStatus
 - Mcs::Usb, 91
- OnTableEntryChanged
 - CWarnerValveControllerDeviceNet, 645
- OnUpdateFirmwareProgress
 - Mcs::Usb, 91
- OnUpdateFirmwareStatusChange
 - Mcs::Usb, 91
- Open
 - Mcs::Usb, 72
- OpenClamp
 - Mcs::Usb, 87
- OpenPipe
 - CGenericDevelopDeviceNet, 157
- OpenPlateClamp
 - CMultiwellDeviceNet, 379
- operator=
 - DeviceIdNet, 664
- OpticalStimulation
 - HeadStageIDType, 677
- Order
 - CCreateFilterNet, 112
 - CFilterPropertyNet, 136
- Output
 - Mcs::Usb, 73
- PacketFrameContextGroup
 - Mcs::Usb, 57, 69, 79, 89
- PatchServAdcModeEnumNet
 - Mcs::Usb, 72
- PatchServer
 - Mcs::Usb, 76
- PathIdent
 - Mcs::Usb, 75
- PatternListStart
 - COctoPotDeviceNet, 389
- PC
 - Mcs::Usb, 90
- PCI
 - Mcs::Usb, 87
- PCX
 - Mcs::Usb, 74
- PeriodicPulse
 - Mcs::Usb, 60
- PeristalticPump
 - Mcs::Usb, 75
- PGA
 - Mcs::Usb, 74
- PIC
 - FirmwareDestinationNames, 676
 - Mcs::Usb, 52
- PIC10
 - Mcs::Usb, 52
- PIC11
 - Mcs::Usb, 52
- PIC12
 - Mcs::Usb, 52
- PIC2
 - FirmwareDestinationNames, 676
 - Mcs::Usb, 52
- PIC3
 - FirmwareDestinationNames, 676
 - Mcs::Usb, 52
- PIC4
 - FirmwareDestinationNames, 676
 - Mcs::Usb, 52
- PIC5
 - Mcs::Usb, 52
- PIC6
 - Mcs::Usb, 52
- PIC7
 - Mcs::Usb, 52
- PIC8
 - Mcs::Usb, 52
- PIC9
 - Mcs::Usb, 52
- PicoAmpere
 - Mcs::Usb, 57
- Plate_24W030MGA
 - Mcs::Usb, 72
- Plate_24W300_30_1152GBA
 - Mcs::Usb, 72
- Plate_24W300_30GBA
 - Mcs::Usb, 72
- Plate_24W300_30GBB
 - Mcs::Usb, 72
- Plate_24W300_30GMA
 - Mcs::Usb, 72
- Plate_24W700_100FMA
 - Mcs::Usb, 72
- Plate_24W700_100FMB
 - Mcs::Usb, 72
- Plate_24W700_100FMC
 - Mcs::Usb, 72
- Plate_24W700_100PBA
 - Mcs::Usb, 72
- Plate_72W500_100FMA
 - Mcs::Usb, 72
- Plate_72W500_100PMA
 - Mcs::Usb, 72

- Plate_96W300_80_1152FMA
 - Mcs::Usb, [72](#)
- Plate_96W400_80_1152FMB
 - Mcs::Usb, [72](#)
- Plate_96W700_100FMA
 - Mcs::Usb, [72](#)
- Plate_96W700_100FMB
 - Mcs::Usb, [72](#)
- Plate_96W700_100GBA
 - Mcs::Usb, [72](#)
- Plate_96W700_100GBB
 - Mcs::Usb, [72](#)
- Plate_96W700_100GBC
 - Mcs::Usb, [72](#)
- Plate_96W700_100GBD
 - Mcs::Usb, [72](#)
- Plate_96W700_100GMA
 - Mcs::Usb, [72](#)
- Plate_Dummy
 - Mcs::Usb, [72](#)
- Plate_Dummy_126
 - Mcs::Usb, [72](#)
- PlateClampEnumNet
 - Mcs::Usb, [72](#)
- PlateClampLockEnumNet
 - Mcs::Usb, [73](#)
- PlusMinus10Volts
 - Mcs::Usb, [51](#)
- PlusMinus2Comma5Volts
 - Mcs::Usb, [51](#)
- PlusMinus5Volts
 - Mcs::Usb, [51](#)
- PollStatusEvent
 - CStimulusFunctionNet, [567](#)
 - CW2100_StimulatorFunctionNet, [611](#)
- PortDirectionEnumNet
 - Mcs::Usb, [73](#)
- Pos900
 - Mcs::Usb, [76](#)
- PositionBase
 - Mcs::Usb, [76](#)
- PositionIIBase
 - Mcs::Usb, [76](#)
- PositionIICentralUnit
 - Mcs::Usb, [76](#)
- PositionImp
 - Mcs::Usb, [76](#)
- PostCommaA
 - CFilterCoefficientsNet::s_FilterAttributesNet, [693](#)
 - Mcs::Usb, [63](#)
- PostCommaB
 - CFilterCoefficientsNet::s_FilterAttributesNet, [693](#)
 - Mcs::Usb, [63](#)
- PowerChip
 - CCMOSMea_FunctionNet, [106](#)
- PowerHS
 - CSCUFunctionNet, [504](#)
- PP_Pump_Mode_Type_EnumNet
 - Mcs::Usb, [73](#)
- PPC
 - Mcs::Usb, [75](#)
- PPCFunction
 - CPPCDeviceNet, [414](#)
- PPS2
 - Mcs::Usb, [75](#)
- PPS4plus1
 - Mcs::Usb, [75](#)
- PPS5
 - Mcs::Usb, [75](#)
- PPS5_DIG
 - Mcs::Usb, [75](#)
- PPS_Function
 - CPPS_DeviceNet, [423](#)
- PreCommaA
 - CFilterCoefficientsNet::s_FilterAttributesNet, [693](#)
 - Mcs::Usb, [63](#)
- PreCommaB
 - CFilterCoefficientsNet::s_FilterAttributesNet, [693](#)
 - Mcs::Usb, [63](#)
- PrepareAndAppendData
 - CStg200xDownloadNet, [552](#)
 - CStimulusFunctionNet, [562](#)
- PrepareAndSendData
 - CStg200xDownloadNet, [553](#)
 - CStimulusFunctionNet, [563](#)
- PrepareChannelData
 - CDigOutStimulatorFunctionNet, [125](#)
- PrepareData
 - CStimulusFunctionNet, [564](#)
 - CW2100_StimulatorFunctionNet, [609](#)
- PrepareDataSync
 - CW2100_StimulatorFunctionNet, [609](#)
- Product
 - CMcsUsbListEntryNet, [304](#)
- ProductIdEnumNet
 - Mcs::Usb, [73](#)
- Program
 - CProgramPressureCurveNet, [430](#)
- PT100
 - Mcs::Usb, [83](#)
- PT1000
 - Mcs::Usb, [83](#)
- PulseGenerator
 - CW2100_FunctionNet, [604](#)
 - Mcs::Usb, [61](#), [70](#), [79](#), [83](#), [89](#)
- PulseGenerator_Mode_EnumNet
 - Mcs::Usb, [76](#)
- PumpOff
 - CRoboFluidDeviceNet, [479](#)
- PumpOn
 - CRoboFluidDeviceNet, [479](#)
- QueryTriggerstatus
 - CStg200xDownloadNet, [554](#)
- raGate
 - Mcs::Usb, [77](#)

raIgnore
 Mcs::Usb, [77](#)
RampStart
 COctoPotDeviceNet, [389](#)
raRestart
 Mcs::Usb, [77](#)
raSingle
 Mcs::Usb, [77](#)
raStop
 Mcs::Usb, [77](#)
rawdata
 Mcs::Usb, [82](#)
RC
 Mcs::Usb, [64](#)
rc100mAh
 Mcs::Usb, [67](#)
rc200mAh
 Mcs::Usb, [67](#)
rc300mAh
 Mcs::Usb, [67](#)
rc30mAh
 Mcs::Usb, [67](#)
rcGreater300mAh
 Mcs::Usb, [67](#)
Read
 CExternDTesterDeviceNet, [128](#)
Read2
 CExternDTesterDeviceNet, [128](#)
ReadBlockFromFlash
 CMcsUsbFactoryNet, [295](#)
ReadBlockFromIFBGlobalEEProm
 CMcsUsbFactoryNet, [295](#)
ReadBlockFromNVMEM
 CMcsUsbFactoryNet, [295](#)
ReadClipping
 CLIH3DeviceNet, [194](#)
ReadEepromRegisterPreconfig
 CMcsUsbNet, [320](#), [321](#)
ReadPipe
 CGenericDevelopDeviceNet, [157](#)
ReadRegister
 CMcsUsbNet, [321](#)
ReadRegister32
 CMcsUsbNet, [321](#)
ReadRegisterTimeSlot
 CMcsUsbNet, [321](#)
ReadUARTData
 CLIH3DeviceNet, [194](#)
Receive
 CSerialPortNet, [508](#)
ReceiveString
 CSerialPortNet, [508](#)
Rectangle
 Mcs::Usb, [84](#)
Ref16
 Mcs::Usb, [77](#)
Ref24
 Mcs::Usb, [77](#)
Ref32
 Mcs::Usb, [77](#)
Ref8
 Mcs::Usb, [77](#)
Reference
 Mcs::Usb, [77](#)
ReferenceElectrodeModeEnumNet
 Mcs::Usb, [76](#)
ReferenceElectrodeSwitchPositionEnumNet
 Mcs::Usb, [77](#)
RegisterHigh
 Mcs::Usb, [60](#)
RegisterLow
 Mcs::Usb, [60](#)
Regular
 Mcs::Usb, [83](#)
RemoveSoftwareKey
 CMcsUsbNet, [321](#)
Renesas
 Mcs::Usb, [86](#)
RescanHeadstage
 CMcsUsbNet, [321](#)
Reserved1
 Mcs::Usb, [83](#)
Reserved2
 Mcs::Usb, [83](#)
Reserved3
 Mcs::Usb, [83](#)
Reserved4
 Mcs::Usb, [83](#)
Reserved5
 Mcs::Usb, [83](#)
ResetAdcOffset
 COctoPotDeviceNet, [389](#)
ResetChannelmap
 CWClassicFunctionNet, [660](#)
ResetDacOffset
 COctoPotDeviceNet, [389](#)
ResetHighpassFilter
 CFilterConfigurationNet, [132](#)
ResetPipe
 CGenericDevelopDeviceNet, [157](#)
ResetStatus
 CStg200xDownloadBasicNet, [546](#)
Retina_LED
 Mcs::Usb, [74](#)
RetriggerActionEnumNet
 Mcs::Usb, [77](#)
RFFunction
 CPositionIIDeviceNet, [410](#)
Rising
 Mcs::Usb, [61](#)
RoboCurrentModeEnumNet
 Mcs::Usb, [77](#)
RoboDacq
 CHiClampDeviceNet, [178](#)
RoboDevice
 CSafeISDeviceNet, [491](#)

- RoboError_AnotherMaster
 - CRoboDeviceNet, [473](#)
- RoboError_Base
 - CRoboDeviceNet, [474](#)
- RoboError_CannotEscapeEndSwitch
 - CRoboDeviceNet, [474](#)
- RoboError_CommandAlreadyInProgress
 - CRoboDeviceNet, [474](#)
- RoboError_CommandNotPossible
 - CRoboDeviceNet, [474](#)
- RoboError_CommunicationTimeout
 - CRoboDeviceNet, [474](#)
- RoboError_DacqNotReady
 - CRoboDeviceNet, [474](#)
- RoboError_DLLMovementTimeout
 - CRoboDeviceNet, [474](#)
- RoboError_FindReferenceMethod
 - CRoboDeviceNet, [474](#)
- RoboError_GilsonCommandPending
 - CRoboDeviceNet, [474](#)
- RoboError_GilsonTimeout
 - CRoboDeviceNet, [474](#)
- RoboError_GilsonWrondID
 - CRoboDeviceNet, [475](#)
- RoboError_McsBus_UnknownCommand
 - CRoboDeviceNet, [475](#)
- RoboError_NoEndSwitch
 - CRoboDeviceNet, [475](#)
- RoboError_NoMoreData
 - CRoboDeviceNet, [475](#)
- RoboError_NoReference
 - CRoboDeviceNet, [475](#)
- RoboError_NoSpeedOrAcceleration
 - CRoboDeviceNet, [475](#)
- RoboError_OverPressure
 - CRoboDeviceNet, [475](#)
- RoboError_ParameterNotAllowed
 - CRoboDeviceNet, [475](#)
- RoboError_PeristalticTimeout
 - CRoboDeviceNet, [475](#)
- RoboError_Phase0OutOfRange
 - CRoboDeviceNet, [475](#)
- RoboError_PollLoopCanceled
 - CRoboDeviceNet, [476](#)
- RoboError_PollLoopCanceledAndStopMovement
 - CRoboDeviceNet, [476](#)
- RoboError_Pressure
 - CRoboDeviceNet, [476](#)
- RoboError_RangeExceeded
 - CRoboDeviceNet, [476](#)
- RoboError_StateChangeNotPossible
 - CRoboDeviceNet, [476](#)
- RoboError_Timeout
 - CRoboDeviceNet, [476](#)
- RoboError_UnknownCommand
 - CRoboDeviceNet, [476](#)
- RoboInject
 - Mcs::Usb, [76](#)
- RoboMainLowLevelCommand
 - CRoboDeviceNet, [477](#)
- RoboMainStatorLowLevelCommand
 - CRoboStatorDeviceNet, [488](#)
- Robooycte2
 - Mcs::Usb, [76](#)
- RoboStatusEvent
 - CRoboDeviceNet, [477](#)
- RoboStatusEventDelegate
 - Mcs::Usb, [91](#)
- Running
 - Mcs::Usb, [81](#)
- RunTable
 - CRoboDacqNet, [455](#)
- s_FilterAttributesNet
 - CFilterCoefficientsNet::s_FilterAttributesNet, [692](#)
- SafeIS
 - Mcs::Usb, [76](#)
- SampleDstSize16
 - Mcs::Usb, [78](#)
- SampleDstSize32
 - Mcs::Usb, [78](#)
- SampleDstSizeNet
 - Mcs::Usb, [78](#)
- SampleRate
 - CCreateFilterNet, [112](#)
- Samplerate
 - CMcsUsbDacqNet, [287](#)
- SampleSize16Signed
 - Mcs::Usb, [78](#)
- SampleSize16Unsigned
 - Mcs::Usb, [78](#)
- SampleSize24Signed
 - Mcs::Usb, [78](#)
- SampleSize24Unsigned
 - Mcs::Usb, [78](#)
- SampleSize32Signed
 - Mcs::Usb, [78](#)
- SampleSize32Unsigned
 - Mcs::Usb, [78](#)
- SampleSize64Signed
 - Mcs::Usb, [78](#)
- SampleSize64Unsigned
 - Mcs::Usb, [78](#)
- SampleSizeNet
 - Mcs::Usb, [78](#)
- SBSVector1
 - Mcs::Usb, [55](#)
- SBSVector2
 - Mcs::Usb, [55](#)
- SBSVector3
 - Mcs::Usb, [55](#)
- SBSVector4
 - Mcs::Usb, [55](#)
- Scale
 - CCreateFilterNet, [112](#)
- ScanForHeadstages
 - CWClassicFunctionNet, [660](#)

- SCU1ElectrodeGroupHS1
 - Mcs::Usb, [79](#)
- SCU1ElectrodeGroupHS2
 - Mcs::Usb, [79](#)
- SCU1ElectrodeGroupHS3
 - Mcs::Usb, [79](#)
- SCU1ElectrodeGroupHS4
 - Mcs::Usb, [79](#)
- SCU2ElectrodeGroupHS1
 - Mcs::Usb, [79](#)
- SCU2ElectrodeGroupHS2
 - Mcs::Usb, [79](#)
- SCU2ElectrodeGroupHS3
 - Mcs::Usb, [79](#)
- SCU2ElectrodeGroupHS4
 - Mcs::Usb, [79](#)
- SCU_HeadstageIdEnumNet
 - Mcs::Usb, [78](#)
- SCUDacqGroupChannelEnumNet
 - Mcs::Usb, [79](#)
- SCUDigitalSourceEnumNet
 - Mcs::Usb, [79](#)
- SelectHeadstage
 - CW2100_FunctionNet, [603](#)
- SelectTimeSlot
 - CW2100_StimulatorFunctionNet, [609](#)
- Send
 - CSerialPortNet, [508](#)
- SendBuffered
 - CGilsonDeviceNet, [164](#)
- SendChannelData
 - CDigOutStimulatorFunctionNet, [125](#)
 - CStg200xDownloadBasicNet, [546](#)
- SendCommand
 - CLIH3DeviceNet, [194](#)
- SendImmediate
 - CGilsonDeviceNet, [164](#)
- SendImmediateGetResponse
 - CGilsonDeviceNet, [164](#)
- SendMultiplexedData
 - CStimulusFunctionNet, [564](#)
- SendPreparedData
 - CStimulusFunctionNet, [564](#)
 - CW2100_StimulatorFunctionNet, [609](#)
- SendSegmentDefine
 - CStg200xDownloadNet, [554](#)
- SendSegmentSelect
 - CStg200xDownloadNet, [554](#)
- SendSegmentStart
 - CStg200xDownloadNet, [555](#)
- SendStart
 - CStg200xBasicNet, [526](#)
 - CStimulusFunctionNet, [564](#)
 - CW2100_StimulatorFunctionNet, [610](#)
- SendStartDacq
 - CMcsUsbDacqNet, [270](#)
- SendStartStgAndDacq
 - CMcsUsbDacqNet, [270](#)
- SendStop
 - CStg200xBasicNet, [526](#)
 - CStimulusFunctionNet, [565](#)
 - CW2100_StimulatorFunctionNet, [610](#)
- SendStopDacq
 - CMcsUsbDacqNet, [271](#)
- SendStopStgAndDacq
 - CMcsUsbDacqNet, [271](#)
- SendStopStgAndDacqWithOptions
 - CMcsUsbDacqNet, [271](#)
- SendSyncData
 - CStg200xDownloadBasicNet, [546](#)
- Sensor
 - CFYIDeviceNet, [144](#)
 - CMeasureTableDeviceNet, [359](#)
 - CPatchServerDeviceNet, [395](#)
- SerialNumber
 - CMcsUsbListEntryNet, [304](#)
 - CMcsUsbNet, [328](#)
- SerialPort
 - CHLADeviceNet, [179](#)
- Set4ADCCatchampAverageShift
 - CMcsBus_SensorNet, [227](#)
- Set4ADCMode
 - CMcsBus_SensorNet, [227](#)
- Set4DAC
 - CMcsBus_SensorNet, [227](#)
- Set_Values
 - CNF_GenDeviceNet, [387](#)
 - CPathIdentDeviceNet, [396](#)
- SetAbsMaxCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [385](#)
- SetAccelerationI
 - CRoboStatorDeviceNet, [486](#)
- SetAccelerationNative
 - CRoboStatorDeviceNet, [486](#)
- SetAccelerationNativeXY
 - CRoboStatorDeviceNet, [486](#)
- SetAccelerationNativeZ
 - CRoboStatorDeviceNet, [486](#)
- SetAccelerationXY
 - CRoboStatorDeviceNet, [486](#)
- SetAccelerationZ
 - CRoboStatorDeviceNet, [486](#)
- SetAccelGyroDesiredRate
 - CW2100_FunctionNet, [603](#)
- SetAccelGyroEnabled
 - CW2100_FunctionNet, [603](#)
- SetAccelRange
 - CW2100_FunctionNet, [603](#)
- SetActiveRunningTableNumber
 - CWarnerValveControllerDeviceNet, [645](#)
- SetADC
 - CWarnerValveControllerDeviceTesterFunctionNet, [656](#)
- SetAdcChannels
 - CSafeISDeviceNet, [490](#)
- SetADCInputOffset

- CCMOSMea_FunctionNet, 106
- SetAdcOffset
 - CLIH3DeviceNet, 194
 - COctoPotDeviceNet, 390
- SetAdcOffsetPermanent
 - CLIH3DeviceNet, 195
- SetAdcSamplePos
 - CSafeISDeviceNet, 490
- SetAirpressureLimit
 - CRoboDeviceNet, 470
- SetAirValve
 - CRoboDeviceNet, 470
- SetAllDigout
 - CRoboDacqNet, 455
- SetAmplificationSwitch
 - COctoPotDeviceNet, 390
- SetAmplitude
 - CChannelTestDeviceNet, 96
- SetAmplitude_nA
 - CTEERFunctionNet, 590
- SetAnalogOutADCRange
 - CSCUFunctionNet, 505
- SetAnalogOutChannel
 - CW2100_FunctionNet, 603
- SetAnalogOutChannels
 - CSCUFunctionNet, 505
- SetAnalogOutDACRange
 - CSCUFunctionNet, 505
- SetAnalogOutFilter
 - CW2100_FunctionNet, 603
- SetAnalogThresholdHigh
 - CWarnerValveControllerDeviceNet, 645
- SetAnalogThresholdLow
 - CWarnerValveControllerDeviceNet, 646
- SetAnalogVoltageRange
 - CPPCFunctionNet, 421
- SetAnalogVoltages
 - CPPS_FunctionNet, 426
- SetAttenuation
 - CChannelTestDeviceNet, 96
- SetAudioChannels
 - CMeaAudioFunctionNet, 333
 - CW2100_FunctionNet, 603
- SetAudioOutDacParameter
 - CLIH3DeviceNet, 195
- SetAutocalibrationDisabled
 - CStg200xBasicNet, 527
- SetAxisConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, 689
- SetAxisLED
 - CRoboocyte2DeviceNet, 482
- SetAxisParametersEeprom
 - CMcsBus_AxisParametersNet, 198, 199
- SetBandwidthByIndex
 - CLintanMea_FunctionNet, 184
- SetBaseFrequency
 - CRFFFunctionNet, 442
- SetBaseSamplerate
 - CCMOSMeaDeviceNet, 110
 - CGrapheneASICDeviceNet, 166
- SetBath
 - CCMOSMea_FunctionNet, 106
- SetBathclamp
 - COctoPotDeviceNet, 390
- SetBathMode
 - CCMOSMea_FunctionNet, 106
- SetBlankingEnable
 - CStg200xBasicNet, 527, 528
- SetBuffer
 - CGenericDevelopDeviceNet, 158
- SetBufferIndex
 - CTEERFunctionNet, 590
- SetBusAddress
 - CMcsBusNet, 238
- SetBusAddressEeprom
 - CMcsBusNet, 238
- SetByteBuffer
 - CGenericDevelopDeviceNet, 158
- SetCalibration
 - CTcxDeviceNet, 580
- SetCardinalDacqSamplerate
 - CInterfaceboardFunctionNet, 187
- SetCardinalStgOutputrate
 - CInterfaceboardFunctionNet, 188
- SetChannel
 - CSw2to64DeviceNet, 568
- SetChannelmap
 - CWClassicFunctionNet, 660
- SetChannels
 - CSw2to64DeviceNet, 569
- SetChannelSwitch
 - COctoPotDeviceNet, 390
- SetChargingMode
 - CMultiBatteryChargerDeviceNet, 369
- SetChargingPCoefficient
 - CMultiBatteryChargerDeviceNet, 370
- SetCheckVoltage
 - CokuvisionStimulatorDeviceNet, 394
- SetClampMode
 - CTEERFunctionNet, 591
 - CWarnerUssingFunctionNet, 626
- SetColorRgb
 - CMultiwellOptoStimFunctionNet, 385
- SetColorStr
 - CMultiwellOptoStimFunctionNet, 385
- SetCommand
 - CMcsBusNet, 239
 - CPedoterDeviceNet, 398
 - CRoboDacqNet, 455
- SetConfiguration
 - CMcsUsbNet, 322
- SetConfigurationBit
 - CRoboDacqNet, 455
- SetConfigurationBitAxc
 - CRoboDacqNet, 456

SetConfigurationBitBlu_Led
CRoboDacqNet, [456](#)

SetConfigurationBitBlu_LedToggleFast
CRoboDacqNet, [456](#)

SetConfigurationBitBlu_LedToggleSlow
CRoboDacqNet, [456](#)

SetConfigurationBitCC_Gen
CRoboDacqNet, [456](#)

SetConfigurationBitCV_Gen
CRoboDacqNet, [456](#)

SetConfigurationBitRC_Gen
CRoboDacqNet, [456](#)

SetConfigurationBitRed_Led
CRoboDacqNet, [456](#)

SetConfigurationBitRed_LedSaturation
CRoboDacqNet, [456](#)

SetConfigurationBitRed_LedToggleFast
CRoboDacqNet, [457](#)

SetConfigurationBitRed_LedToggleSlow
CRoboDacqNet, [457](#)

SetConfigurationBitRelais
CRoboDacqNet, [457](#)

SetConfigurationBitRV_Gen
CRoboDacqNet, [457](#)

SetConfigurationBitStream
CRoboDacqNet, [457](#)

SetConfigurationBitSupply
CRoboDacqNet, [457](#)

SetControllerParams
CTEERFunctionNet, [591](#)

SetCrossTalkOffset
CRoboDacqNet, [457](#)

SetCrossTalkOptimum
CRoboDacqNet, [457](#)

SetCurrentAirvalveLimit
CRoboDeviceNet, [470](#)

SetCurrentAndAir
CRoboDeviceNet, [470](#)

SetCurrentAndAirXY
CRoboStatorDeviceNet, [486](#)

SetCurrentEditTableNumber
CWarnerValveControllerDeviceNet, [646](#)

SetCurrentEnable
CTEERFunctionNet, [591](#)

SetCurrentFactor
COKuvisionStimulatorDeviceNet, [394](#)

SetCurrentMode
CStg200xBasicNet, [528](#)

SetCurrentRangeByIndex
CStg200xBasicNet, [529](#)

SetCycles
CMeaCleanDeviceNet, [336](#)
CMeaCoatDeviceNet, [341](#)

SetD
CTcxDeviceNet, [580](#)

SetDacAmplificationFactor
CStg200xBasicNet, [529](#)

SetDacAutoControl
COctoPotDeviceNet, [390](#)

SetDacIdleValue
CLIH3DeviceNet, [195](#)

SetDacMode
CSafeISDeviceNet, [490](#)

SetDACOffset
CGrapheneFunctionNet, [173](#), [174](#)
COKuvisionStimulatorDeviceNet, [394](#)

SetDacOffset
CDacCalibrationFunctionNet, [114](#)
CLIH3DeviceNet, [195](#)
COctoPotDeviceNet, [390](#)

SetDacOffsetPermanent
CLIH3DeviceNet, [196](#)

SetDacPeriode
CSafeISDeviceNet, [490](#)

SetDacPulseform
CSafeISDeviceNet, [490](#)

SetDacqLegacyMode
CSCUFunctionNet, [505](#)

SetDacRange
CW2100_FunctionNet, [604](#)

SetDACs
CMcsBus_SensorNet, [227](#)

SetDacUseIdleValue
CLIH3DeviceNet, [196](#)

SetDacValue
COctoPotDeviceNet, [390](#)

SetDataMode
CMcsUsbDacqNet, [272](#)

SetDefault
CWarnerValveControllerDeviceNet, [646](#)

SetDestinationSerialNumber
CMcsUsbFactoryNet, [295](#)

SetDetectionThreshold
CMcsBus_SensorNet, [227](#)

SetDevice
CTcxDeviceNet, [580](#)

SetDeviceId
CUsbDeviceConfigurationFunctionNet, [595](#)

SetDeviceList
CPositionImpDeviceNet, [412](#)

SetDeviceName
CUsbDeviceConfigurationFunctionNet, [595](#)

SetDeviceType
CTcxDeviceNet, [580](#)

SetDevname
CTcxDeviceNet, [580](#)

SetDiagnosticMode
CIntanMea_FunctionNet, [184](#)

SetDigitalData
CMeaDigitalDataFunctionNet, [352](#), [353](#)

SetDigitalOut
CMeaDeviceNet, [347](#)

SetDigitalOutPortInvert
CWarnerValveControllerDeviceNet, [646](#)

SetDigitalOutPortValve
CWarnerValveControllerDeviceNet, [647](#)

- SetDigitalPortDirection
 - CWarnerValveControllerDeviceNet, [647](#)
- SetDigitalSource
 - CMcsUsbDacqNet, [272–274](#)
- SetDigitalStimulatorTrigger
 - CW2100_StimulatorFunctionNet, [610](#)
- SetDigitalStimulatorTriggerSlope
 - CW2100_StimulatorFunctionNet, [610](#)
- SetDigout
 - CFluidControlDeviceNet, [141](#)
 - CRoboDacqNet, [457](#)
- SetDigoutMode
 - CStg200xBasicNet, [529](#)
- SetDigOutState
 - CLIH3DeviceNet, [196](#)
- SetDigoutValue
 - CStg200xBasicNet, [529](#)
- SetDIO
 - CMcsBus_FYIExtensionNet, [201](#)
- SetDischargeCurrentSetPoint
 - CMultiBatteryChargerDeviceNet, [370](#)
- SetDisplayMode
 - CWarnerValveControllerDeviceNet, [647](#)
- SetDisplayText
 - CRoboDacqNet, [458](#)
- SetDownsampleFactor
 - CRoboDacqNet, [458](#)
- SetDSPHighPassByIndex
 - CIntanMea_FunctionNet, [184](#)
- SetDuration
 - CMeaCoatDeviceNet, [342](#)
- SetEepromPage
 - CLIH3DeviceNet, [196](#)
- SetElectrodeDacMux
 - CStg200xBasicNet, [530](#), [531](#)
- SetElectrodeEnable
 - CStg200xBasicNet, [532](#), [533](#)
- SetElectrodeMode
 - CStg200xBasicNet, [534](#), [535](#)
- SetEnableAmplifierProtectionSwitch
 - CStg200xBasicNet, [536](#), [537](#)
- SetEnableHeaterLimit
 - CTcxDeviceNet, [580](#)
- SetEnablePulse
 - CWarnerUssingFunctionNet, [626](#)
- SetEnableThermocouple
 - CTcxDeviceNet, [581](#)
- SetExternalElectrodeEnable
 - CStg200xBasicNet, [537](#), [538](#)
- SetExternalLED
 - CTEERFunctionNet, [591](#)
- SetFAmplification
 - CStg200xBasicNet, [538](#)
- SetFilter
 - CRoboDacqNet, [458](#)
- SetFilterCoeffs
 - CRoboDacqNet, [458](#)
- SetFilterParameter
 - CFilterConfigurationNet, [132](#)
 - CFilterConfigurationRegisterNet, [134](#)
- SetFilterParameterPermanent
 - CFilterConfigurationNet, [132](#)
 - CFilterConfigurationRegisterNet, [134](#)
- SetFilterParametersHeadstage
 - CWClassicFunctionNet, [660](#)
- SetFinalDischargeVoltage
 - CMultiBatteryChargerDeviceNet, [370](#)
- SetFrequency
 - CChannelTestDeviceNet, [96](#)
 - CRadioControlledDevicesNet, [435](#)
- SetGain
 - CPgaDeviceNet, [401](#)
- SetGate
 - CCMOSMea_FunctionNet, [106](#)
- SetGateFloating
 - CCMOSMea_FunctionNet, [106](#)
- SetGateToVOP
 - CCMOSMea_FunctionNet, [106](#)
- SetGlobalRepeat
 - CDigOutStimulatorFunctionNet, [125](#)
- SetGyroRange
 - CW2100_FunctionNet, [604](#)
- SetHasChecksum
 - CWClassicFunctionNet, [661](#)
- SetHeadstage
 - CStg200xBasicNet, [538](#)
- SetHeadstageOnOff
 - CW2100_FunctionNet, [604](#)
 - CWClassicFunctionNet, [661](#)
- SetHeadstagePowerStateAtStart
 - CSCUFunctionNet, [506](#)
- SetHeadstageSamplingActive
 - CW2100_FunctionNet, [604](#)
- SetHeadstageToSleep
 - CW2100_FunctionNet, [604](#)
- SetHeaterLimit
 - CTcxDeviceNet, [581](#)
- SetHighCurrentMode
 - CWarnerUssingFunctionNet, [626](#)
- SetHighpassFilterEnable
 - CFilterConfigurationNet, [133](#)
- SetHWConfig
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetHWRevision
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetHWRevisionEeprom
 - CMcsBusNet, [239](#)
- SetHWSelectedChannels
 - CWClassicFunctionNet, [661](#)
- SetI
 - CTcxDeviceNet, [581](#)
- SetI Clamp
 - CRoboDacqNet, [458](#)
- SetICoeff
 -

- CRobo_FYITemp_FunctionNet, [446](#)
- SetICOffset
 - CRoboDacqNet, [458](#)
- SetIdleModeOffset
 - CWarnerUssingFunctionNet, [627](#)
- SetIGain
 - CRoboDacqNet, [458](#)
- SetImpedanceTestFrequency
 - CMealImpedanceDeviceNet, [358](#)
- SetImpId
 - CPositionImpDeviceNet, [412](#)
- SetImplantCurrentSetpoint
 - CPositionIIDeviceNet, [407](#)
- SetInMovement
 - CRoboDeviceNet, [470](#)
- SetIntanRegister
 - CIntanMea_FunctionNet, [184](#)
- SetIntBuffer
 - CGenericDevelopDeviceNet, [158](#)
- SetIO
 - CWarnerValveControllerDeviceTesterFunctionNet, [656](#)
- SetIODirection
 - CWarnerValveControllerDeviceTesterFunctionNet, [657](#)
- SetIoVoltage
 - CInterfaceboard2FunctionNet, [186](#)
- SetLatency
 - CMcsBus_SensorNet, [227](#)
- SetLayoutConfiguration
 - CMEA2100x256FunctionNet, [330](#)
- SetLED
 - CRetinalLedDeviceNet, [437](#)
- SetLEDSwitch
 - CMcsBus_ExtensionNet, [200](#)
- SetLength
 - CRobo_FYIProgram_FunctionNet, [444](#)
- SetLiquidResistance
 - CTEERFunctionNet, [592](#)
 - CWarnerUssingFunctionNet, [627](#)
- SetListmodeIndexRange
 - CStg200xBasicNet, [538](#)
- SetListmodeTriggerSource
 - CStg200xBasicNet, [539](#)
- SetLowCurrentMode
 - CWarnerUssingFunctionNet, [627](#)
- SetLumi
 - CRetinalLedDeviceNet, [437](#)
- SetMaxCurrent
 - CMeaCoatDeviceNet, [342](#)
- SetMaxDurationHighCurrentInMicroSec
 - CMultiwellOptoStimFunctionNet, [385](#)
- SetMaxDutyCycleHighCurrent
 - CMultiwellOptoStimFunctionNet, [386](#)
- SetMaxHeaterPowerMultiwell
 - CTcxDeviceNet, [581](#)
- SetMaxNoPressure
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMaxNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMaxP
 - CTcxDeviceNet, [581](#)
- SetMaxPower
 - CokuvisionStimulatorDeviceNet, [394](#)
 - CRobo_FYITemp_FunctionNet, [446](#)
- SetMaxPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMaxVoltage
 - CMeaCleanDeviceNet, [337](#)
 - CokuvisionStimulatorDeviceNet, [394](#)
- SetMCAcceleration
 - CMcsBus_MotorControlNet, [211](#)
- SetMCAccelerationEeprom
 - CMcsBus_MotorControlNet, [212](#)
- SetMCAccelerationShortCommand
 - CMcsBus_MotorControlNet, [212](#)
- SetMCAxisRevisionEeprom
 - CMcsBus_MotorControlNet, [212](#)
- SetMCBreakCurrent
 - CMcsBus_MotorControlNet, [212](#)
- SetMCBreakCurrentEeprom
 - CMcsBus_MotorControlNet, [212](#)
- SetMCConfig
 - CMcsBus_MotorControlNet, [212](#)
- SetMCConfigEeprom
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrent
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrentEeprom
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrentMode
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrentModeEeprom
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrentModeShortCommand
 - CMcsBus_MotorControlNet, [213](#)
- SetMCCurrentPosition
 - CMcsBus_MotorControlNet, [214](#)
- SetMCCurrentShortCommand
 - CMcsBus_MotorControlNet, [214](#)
- SetMCMaxAcceleration
 - CMcsBus_MotorControlNet, [214](#)
- SetMCMaxAccelerationEeprom
 - CMcsBus_MotorControlNet, [214](#)
- SetMCMaxCurrent
 - CMcsBus_MotorControlNet, [214](#)
- SetMCMaxCurrentEeprom
 - CMcsBus_MotorControlNet, [214](#)
- SetMCMaxSpeed
 - CMcsBus_MotorControlNet, [215](#)
- SetMCMaxSpeedEeprom
 - CMcsBus_MotorControlNet, [215](#)

- SetMCMaxTravel
 - CMcsBus_MotorControlNet, [215](#)
- SetMCMaxTravelEeprom
 - CMcsBus_MotorControlNet, [215](#)
- SetMCMaxTravelShortCommand
 - CMcsBus_MotorControlNet, [215](#)
- SetMCNewPosition
 - CMcsBus_MotorControlNet, [215](#)
- SetMCOutputOnOff
 - CMcsBus_MotorControlNet, [216](#)
- SetMCReference
 - CMcsBus_MotorControlNet, [216](#)
- SetMCReferenceCurrent
 - CMcsBus_MotorControlNet, [216](#)
- SetMCReferenceCurrentEeprom
 - CMcsBus_MotorControlNet, [216](#)
- SetMCRegulatorGain
 - CMcsBus_MotorControlNet, [216](#)
- SetMCRegulatorGainEeprom
 - CMcsBus_MotorControlNet, [216](#)
- SetMCRotation
 - CMcsBus_MotorControlNet, [217](#)
- SetMCScalingFactor
 - CMcsBus_MotorControlNet, [217](#)
- SetMCScalingFactorEeprom
 - CMcsBus_MotorControlNet, [217](#)
- SetMCSpeed
 - CMcsBus_MotorControlNet, [217](#)
- SetMCSpeedEeprom
 - CMcsBus_MotorControlNet, [217](#)
- SetMCSpeedShortCommand
 - CMcsBus_MotorControlNet, [217](#)
- SetMCSpeedUnitEeprom
 - CMcsBus_MotorControlNet, [218](#)
- SetMCStandbyCurrent
 - CMcsBus_MotorControlNet, [218](#)
- SetMCStandbyCurrentEeprom
 - CMcsBus_MotorControlNet, [218](#)
- SetMCStandbyTime
 - CMcsBus_MotorControlNet, [218](#)
- SetMCStandbyTimeEeprom
 - CMcsBus_MotorControlNet, [218](#)
- SetMeasurementMode
 - CStg200xBasicNet, [540](#)
- SetMinimalThreshold
 - CMcsBus_SensorNet, [227](#)
- SetMinNoPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMinPressure
 - CRoboDeviceNet, [471](#)
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMinPressureWaitTime
 - CRoboDeviceNet::RoboMainLowLevelCommands, [689](#)
- SetMinVoltage
 - CMeaCleanDeviceNet, [337](#)
- SetModeSelect
 - CPulseGeneratorFunctionNet, [433](#)
- SetMovePump
 - CMcsBus_SensorNet, [228](#)
- SetMultiHeadstageMode
 - CW2100_FunctionNet, [604](#)
- SetNanoVoltsPerKelvin
 - CMcsBus_TempSensorNet, [231](#)
- SetNeurochipMemoryData
 - CCMOSMea_FunctionNet, [106](#)
- SetNoFilterCoeffs
 - CRoboDacqNet, [458](#)
- SetNumberOfAnalogChannels
 - CMeaDeviceNet, [348](#)
- SetNumberOfChannels
 - CMeaDeviceNet, [349](#)
 - COctoPotDeviceNet, [390](#)
- SetOffsetCurrent
 - CMeaCoatDeviceNet, [342](#)
- SetOnOff
 - CTcxDeviceNet, [581](#)
- SetOutputMap
 - CStg200xDownloadNet, [555](#)
- SetOutputRate
 - COctoPotDeviceNet, [390](#)
 - CStg200xBasicNet, [540](#)
- SetP
 - CTcxDeviceNet, [582](#)
- SetParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetPattern
 - CMeaSwitchDeviceNet, [361](#)
- SetPatternBool
 - CMeaSwitchDeviceNet, [361](#)
- SetPatternListEntry
 - COctoPotDeviceNet, [391](#)
- SetPauseDuration
 - CMeaCoatDeviceNet, [342](#)
- SetPCoeff
 - CRobo_FYITemp_FunctionNet, [446](#)
- SetPeriod
 - CPulseGeneratorFunctionNet, [433](#)
- SetPeriod_us
 - CTEERFunctionNet, [592](#)
- SetPermanentCurrentInMicroAmp
 - CMultiwellOptoStimFunctionNet, [386](#)
- SetPersistence
 - CRetinaLedDeviceNet, [437](#)
- SetPGain
 - CRoboDacqNet, [458](#)
- SetPidParameter
 - COctoPotDeviceNet, [391](#)
- SetPiezoState
 - CMcsBus_SensorNet, [228](#)
- SetPlateMux
 - CMultiwellDeviceNet, [379](#), [380](#)
- SetPlateType

- CMultiwellDeviceNet, [380](#)
- SetPoti
 - CMcsUsbDacqNet, [274](#)
- SetPowerMuxPlate
 - CMultiwellDeviceNet, [380](#)
- SetPowerStrength
 - CPositionIIDeviceNet, [407](#)
- SetPressureOffset
 - CMcsBus_SensorNet, [228](#)
 - CPPCFunctionNet, [421](#)
- SetPressureRange
 - CPPCFunctionNet, [421](#)
- SetPulse
 - CWarnerUssingFunctionNet, [627](#)
- SetPulseform
 - CokuvisionStimulatorDeviceNet, [394](#)
- SetPulseLength
 - CPulseGeneratorFunctionNet, [433](#)
- SetPumpCouple
 - CPPS_FunctionNet, [427](#)
- SetPumpEnableSpeedRatio
 - CPPS_FunctionNet, [427](#)
- SetPumpFastOnOff
 - CPPS_FunctionNet, [427](#)
- SetPumpFastSpeed
 - CPPS_FunctionNet, [427](#)
- SetPumpFunctionSpeeds
 - CPPS_FunctionNet, [427](#)
- SetPumpManualOnOff
 - CPPS_FunctionNet, [427](#)
- SetPumpMaxSpeed
 - CPPS_FunctionNet, [427](#)
- SetPumpModeType
 - CPPCFunctionNet, [422](#)
 - CPPS_FunctionNet, [427](#)
- SetPumpSpeed
 - CRoboFluidDeviceNet, [479](#)
- SetPumpSpeedRatio
 - CPPS_FunctionNet, [428](#)
- SetPumpSpeedUnit
 - CPPCFunctionNet, [422](#)
 - CPPS_FunctionNet, [428](#)
- SetPWM
 - CFluidControlDeviceNet, [141](#)
- SetRampParameter
 - COctoPotDeviceNet, [391](#)
- SetRatedCapacity
 - CMultiBatteryChargerDeviceNet, [370](#)
- SetRatedCapacityVolatile
 - CMultiBatteryChargerDeviceNet, [372](#)
- SetRecordingNumber
 - CRoboDacqNet, [459](#)
- SetReferenceElectrodeMode
 - CSCUFunctionNet, [506](#)
- SetReferenceElectrodeSwitchState
 - CSCUFunctionNet, [506](#)
- SetRegionOfInterests
 - CCMOSMeaDeviceNet, [110](#)
- CGrapheneASICDeviceNet, [166](#)
- SetRegulationTimeouts
 - CMcsBus_SensorNet, [228](#)
- SetRegulatorFactor
 - CMcsBus_SensorNet, [228](#)
- SetRegulatorOnOff
 - CMcsBus_SensorNet, [228](#)
 - CRobo_FYITemp_FunctionNet, [446](#)
- SetRepeat
 - CRetinaLedDeviceNet, [438](#)
- SetRepeats
 - CProgramPressureCurveNet, [430](#)
- SetResetFilter
 - CWClassicFunctionNet, [661](#)
- SetRFFrequency
 - CPositionImpDeviceNet, [413](#)
- SetRFFrequencyHeadstage
 - CWClassicFunctionNet, [661](#)
- SetRFFrequencyReceiver
 - CWClassicFunctionNet, [661](#)
- SetRFFrequencyReceiverEeprom
 - CWClassicFunctionNet, [661](#)
- SetRFLostBehaviour
 - CWClassicFunctionNet, [661](#)
- SetRFPower
 - CWClassicFunctionNet, [662](#)
- SetRotatePump
 - CMcsBus_SensorNet, [228](#)
- SetRTC
 - CokuvisionStimulatorDeviceNet, [394](#)
 - CPositionIIDeviceNet, [407](#)
- SetSampleInterval
 - CLIH3DeviceNet, [197](#)
- SetSamplePeriode
 - CMcsBus_SensorNet, [229](#)
- SetSamplerate
 - CMcsUsbDacqNet, [274](#)
- SetScreen
 - CRoboDacqNet, [459](#)
- SetSearchReferenceFastAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetSearchReferenceFastSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetSearchReferenceFineAccel
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetSearchReferenceFineSpeed
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetSearchReferenceMethod
 - CRoboDeviceNet::RoboMainLowLevelCommands, [690](#)
- SetSearchReferenceMoveOut
 - CRoboDeviceNet::RoboMainLowLevelCommands, [691](#)
- SetSearchReferenceOffsetPos

- CRoboDeviceNet::RoboMainLowLevelCommands, 691
- SetSelectedChannels
 - CMcsUsbDacqNet, 275–277
 - CW2100_FunctionNet, 604
- SetSelectedChannelsQueue
 - CMcsUsbDacqNet, 277–279
- SetSelectedData
 - CMcsUsbDacqNet, 279–281
- SetSelectedHeadstage
 - CWClassicFunctionNet, 662
- SetSensorType
 - CTcxDeviceNet, 582
- SetSerialNumberHeadstage
 - CWClassicFunctionNet, 662
- SetSetpoint
 - CTcxDeviceNet, 582
- SetShortBuffer
 - CGenericDevelopDeviceNet, 159
- SetSimulation
 - CRoboDacqNet, 459
- SetSineParameter
 - COctoPotDeviceNet, 391
- SetSingleHeater
 - CMcsBus_FYIExtensionNet, 201
- SetSingleValve
 - CFluidControlDeviceNet, 141
 - CRoboFluidDeviceNet, 479
- SetSlope
 - CMeaCleanDeviceNet, 337
 - CMeaCoatDeviceNet, 343
- SetSoftwareKey
 - CMcsUsbNet, 322
- SetSollPressure
 - CMcsBus_SensorNet, 229
- SetSollTemp
 - CRobo_FYITemp_FunctionNet, 446
- SetSourceBulk
 - CCMOSMea_FunctionNet, 106
- SetSourceDrain
 - CCMOSMea_FunctionNet, 107
- SetSourceGate
 - CCMOSMea_FunctionNet, 107
- SetSpeedI
 - CRoboStatorDeviceNet, 487
- SetSpeedNativeI
 - CRoboStatorDeviceNet, 487
- SetSpeedNativeXY
 - CRoboStatorDeviceNet, 487
- SetSpeedNativeZ
 - CRoboStatorDeviceNet, 487
- SetSpeedXY
 - CRoboStatorDeviceNet, 487
- SetSpeedZ
 - CRoboStatorDeviceNet, 487
- SetStartTriggerSlope
 - CDigOutStimulatorFunctionNet, 126
- SetStateDebugData
 - CPositionIIDeviceNet, 409
- SetStateEventData
 - CPositionIIDeviceNet, 409
- SetStgProgramInfo
 - CStg200xBasicNet, 540
- SetStimulusSites
 - CCMOSMea_FunctionNet, 107
- SetStopTriggerSlope
 - CDigOutStimulatorFunctionNet, 126
- SetStringFormat
 - CMcsUsbListEntryNet, 303
 - CMcsUsbListNet, 307
- SetSubChannel
 - CMcsBus_MotorControlNet, 218
- SetSwitches
 - CSafeISDeviceNet, 491
- SetSyncoutMap
 - CStg200xBasicNet, 541
- SetTableName
 - CWarnerValveControllerDeviceNet, 647
- SetTablepointer
 - CRetinaLedDeviceNet, 438
- SetTableStep
 - CWarnerValveControllerDeviceNet, 648
- SetTableStepAll
 - CWarnerValveControllerDeviceNet, 648
- SetTestMode
 - CRFFFunctionNet, 442
- SetThermocoupleNanovoltPerKelvin
 - CFluidControlDeviceNet, 143
 - CTcxDeviceNet, 582
- SetThermoOffset
 - CMcsBus_TempSensorNet, 231
- Settings
 - Mcs::Usb, 90
- SetTouchPadEnable
 - CMultiwellDeviceNet, 381
- SetTransformer
 - CMeFunctionNet, 364
- SetTrigger
 - CRetinaLedDeviceNet, 438
 - CWarnerValveControllerDeviceTesterFunctionNet, 657
- SetTriggerMaskValue
 - CMeaDeviceNet, 349
 - CRoboDacqNet, 459
- SetTriggerPeriod
 - CMeaDeviceNet, 350
- SetTriggerSource
 - CStg200xBasicNet, 541
- SetTriggerSyncDirection
 - CWarnerValveControllerDeviceTesterFunctionNet, 657
- SetUByteBuffer
 - CGenericDevelopDeviceNet, 160
- SetUClamp
 - CRoboDacqNet, 459
- SetUCOffset

- CRoboDacqNet, 459
- SetUIntBuffer
 - CGenericDevelopDeviceNet, 160
- SetupGroupDacqQueue
 - CMcsUsbDacqNet, 281
- SetupRetriggerMode
 - CStg200xDownloadBasicNet, 547
- SetupTrigger
 - CStg200xDownloadBasicNet, 548
 - CStimulusFunctionNet, 565
- SetupTriggerSingle
 - CStg200xDownloadBasicNet, 548
 - CStimulusFunctionNet, 566
- SetUseBubble
 - CPPS_FunctionNet, 428
- SetUserParameter
 - CRoboDeviceNet::RoboMainLowLevelCommands, 691
- SetUShortBuffer
 - CGenericDevelopDeviceNet, 161
- SetUOffset
 - CRoboDacqNet, 459
- SetValue
 - CGenericDevelopDeviceNet, 161
- SetValve
 - CFluidControlDeviceNet, 143
 - CRoboFluidDeviceNet, 479
- SetValve1
 - CRobo_FYIProgram_FunctionNet, 444
- SetValve2
 - CRobo_FYIProgram_FunctionNet, 444
- SetValveActive
 - CPPCFunctionNet, 422
 - CWarnerValveControllerDeviceNet, 648
- SetValveCurrent
 - CWarnerValveControllerDeviceNet, 648
- SetValveDigitalInInvert
 - CWarnerValveControllerDeviceNet, 649
- SetValveDigitalInPort
 - CWarnerValveControllerDeviceNet, 649
- SetValveLedOn
 - CWarnerValveControllerDeviceNet, 649
- SetValveManualGroup
 - CWarnerValveControllerDeviceNet, 649
- SetValveManualState
 - CWarnerValveControllerDeviceNet, 650
- SetValveMode
 - CWarnerValveControllerDeviceNet, 650
- SetValves
 - CMcsBus_FYIExtensionNet, 201
- SetValvesActiveMap
 - CWarnerValveControllerDeviceNet, 650
- SetValvesManualStateMap
 - CWarnerValveControllerDeviceNet, 650
- SetValveTableEntry
 - CWarnerValveControllerDeviceNet, 650
- SetVdsVgs
 - CGrapheneFunctionNet, 174
- SetVdVs
 - CGrapheneFunctionNet, 175
- SetVdVsDAC
 - CGrapheneFunctionNet, 175
- SetVelocityI
 - CRoboStatorDeviceNet, 487
- SetVelocityXY
 - CRoboStatorDeviceNet, 487
- SetVelocityZ
 - CRoboStatorDeviceNet, 487
- SetVMMaxNegativeCurrent
 - CMcsBus_VoltageModeNet, 234
- SetVMMaxNegativeCurrentEeprom
 - CMcsBus_VoltageModeNet, 234
- SetVMMaxNegativeVoltage
 - CMcsBus_VoltageModeNet, 234
- SetVMMaxNegativeVoltageEeprom
 - CMcsBus_VoltageModeNet, 234
- SetVMMaxPositiveCurrent
 - CMcsBus_VoltageModeNet, 235
- SetVMMaxPositiveCurrentEeprom
 - CMcsBus_VoltageModeNet, 235
- SetVMMaxPositiveVoltage
 - CMcsBus_VoltageModeNet, 235
- SetVMMaxPositiveVoltageEeprom
 - CMcsBus_VoltageModeNet, 235
- SetVMOOutputOnOff
 - CMcsBus_VoltageModeNet, 235
- SetVMVoltage
 - CMcsBus_VoltageModeNet, 235
- SetVolatileClampOffset
 - CMultiwellDeviceNet, 381
- SetVoltage12VLimit
 - CRoboDeviceNet, 471
- SetVoltage5VLimit
 - CRoboDeviceNet, 471
- SetVoltageAirvalveLimit
 - CRoboDeviceNet, 471
- SetVoltageClampControllerParam_D
 - CWarnerUssingFunctionNet, 628
- SetVoltageClampControllerParam_I
 - CWarnerUssingFunctionNet, 628
- SetVoltageClampControllerParam_P
 - CWarnerUssingFunctionNet, 628
- SetVoltageMode
 - CStg200xBasicNet, 541
- SetVoltageRange
 - CGrapheneFunctionNet, 176
- SetVoltageRangeByIndex
 - CMcsUsbDacqNet, 281
- SetVoltageRangeInMicroVolt
 - CMcsUsbDacqNet, 282
- SetVoltageResolution
 - CGrapheneFunctionNet, 176
- SetVoltageRs485ALimit
 - CRoboDeviceNet, 471
- SetVoltageRs485BLimit
 - CRoboDeviceNet, 471

- SetVoltageValvesLimit
 - CRoboDeviceNet, [471](#)
- SetWaveform
 - CChannelTestDeviceNet, [96](#)
 - CTEERFunctionNet, [592](#)
- SetWaveLengthInNanometer
 - CMultiwellOptoStimFunctionNet, [386](#)
- SetWorkingFrequency
 - CRFFFunctionNet, [442](#)
- SetWPADebugMode
 - CWClassicFunctionNet, [662](#)
- SetWPAType
 - CWClassicFunctionNet, [662](#)
- SetXGain
 - CRoboDacqNet, [459](#)
- Sideband
 - CStimulusFunctionNet::SidebandData, [694](#)
- SidebandData
 - CStimulusFunctionNet::SidebandData, [694](#)
- Signed_16bit
 - Mcs::Usb, [58](#)
- Signed_24bit
 - Mcs::Usb, [58](#)
- Signed_32bit
 - Mcs::Usb, [58](#)
- Sine
 - Mcs::Usb, [84](#)
- SineStart
 - COctoPotDeviceNet, [391](#)
- SingleWell
 - Mcs::Usb, [82](#)
- SixWell
 - Mcs::Usb, [82](#)
- size
 - DigitalSource< digitalsourceenum >, [665](#)
 - DigitalSourceGeneral, [667](#)
- SmartImplant
 - Mcs::Usb, [76](#)
- SN
 - HeadStageIDType, [679](#)
- SOFAndCTRLword
 - Mcs::Usb, [56](#)
- Software
 - Mcs::Usb, [64](#)
- SoftwareDongle
 - Mcs::Usb, [75](#)
- Source
 - DigitalSource< digitalsourceenum >, [666](#)
 - DigitalSourceGeneral, [667](#)
- Standby
 - Mcs::Usb, [77](#), [87](#)
- Start
 - CMeaCleanDeviceNet, [337](#)
 - CMeaCoatDeviceNet, [343](#)
 - CRobo_FYIProgram_FunctionNet, [444](#)
 - Mcs::Usb, [61](#)
- StartDacq
 - CMcsUsbDacqNet, [282](#), [283](#)
- StartInternalCalibration
 - CTEERFunctionNet, [592](#)
- StartLoop
 - CMcsUsbDacqNet, [284](#), [285](#)
- StartMCMovement
 - CMcsBus_MotorControlNet, [219](#)
- StartMeasurement
 - CMeaImpedanceDeviceNet, [358](#)
- StartPoll
 - CStimulusFunctionNet, [566](#)
 - CW2100_StimulatorFunctionNet, [610](#)
- StartQueue
 - CRoboDeviceNet, [471](#)
- StartSampling
 - CTEERFunctionNet, [592](#)
- StartSync
 - CMcsBus_SensorNet, [229](#)
- State
 - HeadStageIDTypeState, [681](#)
 - Mcs::Usb, [87](#)
- Status
 - CUsbExceptionNet, [596](#)
- Status_AlreadyConfigured
 - CMcsUsbNet, [324](#)
- Status_BadStartFrame
 - CMcsUsbNet, [324](#)
- Status_Btstuff
 - CMcsUsbNet, [324](#)
- Status_BufferOverrun
 - CMcsUsbNet, [324](#)
- Status_BufferUnderrun
 - CMcsUsbNet, [324](#)
- Status_Canceled
 - CMcsUsbNet, [324](#)
- Status_Canceling
 - CMcsUsbNet, [325](#)
- Status_ConnectedPipes
 - CMcsUsbNet, [325](#)
- Status_ControlNotOwned
 - CMcsUsbNet, [325](#)
- Status_Crc
 - CMcsUsbNet, [325](#)
- Status_DataOverrun
 - CMcsUsbNet, [325](#)
- Status_DataToggleMismatch
 - CMcsUsbNet, [325](#)
- Status_DataUnderrun
 - CMcsUsbNet, [325](#)
- Status_DeviceLocked
 - CMcsUsbNet, [325](#)
- Status_DeviceNotFound
 - CMcsUsbNet, [325](#)
- Status_DeviceRemoved
 - CMcsUsbNet, [325](#)
- Status_DevNotResponding
 - CMcsUsbNet, [325](#)
- Status_EndpointHalted
 - CMcsUsbNet, [326](#)

- Status_ErrorBusy
 - CMcsUsbNet, [326](#)
- Status_ErrorShortTransfer
 - CMcsUsbNet, [326](#)
- Status_Fifo
 - CMcsUsbNet, [326](#)
- Status_FrameControlOwned
 - CMcsUsbNet, [326](#)
- Status_InternalHcError
 - CMcsUsbNet, [326](#)
- Status_InvalidDeviceHandle
 - CMcsUsbNet, [326](#)
- Status_InvalidHandle
 - CMcsUsbNet, [326](#)
- Status_InvalidParameter
 - CMcsUsbNet, [326](#)
- Status_InvalidPipeHandle
 - CMcsUsbNet, [326](#)
- Status_InvalidUrbFunction
 - CMcsUsbNet, [327](#)
- Status_IoPending
 - CMcsUsbNet, [327](#)
- Status_IoTimeout
 - CMcsUsbNet, [327](#)
- Status_IsochRequestFailed
 - CMcsUsbNet, [327](#)
- Status_LastUsbErrorMismatch
 - CMcsUsbNet, [327](#)
- Status_NoBandwidth
 - CMcsUsbNet, [327](#)
- Status_NoMemory
 - CMcsUsbNet, [327](#)
- Status_NoSuchDevice
 - CMcsUsbNet, [327](#)
- Status_NotAccessed
 - CMcsUsbNet, [327](#)
- Status_NotSupported
 - CMcsUsbNet, [327](#)
- Status_PidCheckFailure
 - CMcsUsbNet, [327](#)
- Status_PipeNotLinked
 - CMcsUsbNet, [328](#)
- Status_RequestFailed
 - CMcsUsbNet, [328](#)
- Status_RequestMutexFailed
 - CMcsUsbNet, [328](#)
- Status_RequestMutexTimeout
 - CMcsUsbNet, [328](#)
- Status_Stall
 - CMcsUsbNet, [328](#)
- Status_Unconfigured
 - CMcsUsbNet, [328](#)
- Status_UnexpectedPid
 - CMcsUsbNet, [328](#)
- STG
 - Mcs::Usb, [74](#)
- Stg1
 - Mcs::Usb, [62](#)
- STG1DACSignalGroup
 - Mcs::Usb, [57](#), [69](#), [79](#)
- STG1SidebandsGroup
 - Mcs::Usb, [57](#), [69](#), [79](#)
- STG1TriggerStatusGroup
 - Mcs::Usb, [57](#), [69](#), [79](#)
- Stg2
 - Mcs::Usb, [62](#)
- Stg200xDigoutModeEnumNet
 - Mcs::Usb, [81](#)
- Stg200xPollStatusEvent
 - CStg200xDownloadNet, [556](#)
- Stg200xSegmentFlagsEnumNet
 - Mcs::Usb, [81](#)
- Stg200xTriggerStatusEnumNet
 - Mcs::Usb, [81](#)
- STG2DACSignalGroup
 - Mcs::Usb, [69](#), [79](#)
- STG2SidebandsGroup
 - Mcs::Usb, [69](#), [79](#)
- STG2TriggerStatusGroup
 - Mcs::Usb, [69](#), [79](#)
- Stg3
 - Mcs::Usb, [62](#)
- STG3008_FA
 - Mcs::Usb, [74](#)
- STG4002
 - Mcs::Usb, [74](#)
- STG4002_opto
 - Mcs::Usb, [74](#)
- STG4004
 - Mcs::Usb, [74](#)
- STG4004_opto
 - Mcs::Usb, [74](#)
- STG4008
 - Mcs::Usb, [74](#)
- STG4008_opto
 - Mcs::Usb, [74](#)
- STG400x
 - Mcs::Usb, [74](#)
- STG400x_opto
 - Mcs::Usb, [74](#)
- STG5
 - Mcs::Usb, [74](#)
- STG_DestinationEnumNet
 - Mcs::Usb, [82](#)
- StgListModeTrigger
 - Mcs::Usb, [62](#)
- StgStatusNet, [695](#)
 - FromIntPtr, [695](#)
 - FromPtr, [695](#)
 - ListOfChangedTriggers, [695](#)
 - TiggerStatus, [695](#)
- StgTrigger
 - Mcs::Usb, [62](#)
- StillConnected
 - CRadioControlledDevicesNet, [435](#)
- Stimulation

- Mcs::Usb, [54](#)
- StimulationLayoutConfigurationEnumNet
 - Mcs::Usb, [82](#)
- Stimulator
 - CW2100_FunctionNet, [605](#)
- Stimulus
 - CCMOSMeaDeviceNet, [110](#)
 - CStg200xDownloadBasicNet, [549](#)
- StimulusDeviceDataAndUnrolledData
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- StimulusFunction
 - CLIH3DeviceNet, [197](#)
- StimulusParameters
 - HeadStageIDType, [679](#)
- Stop
 - CMeaCleanDeviceNet, [337](#)
 - CMeaCoatDeviceNet, [343](#)
 - Mcs::Usb, [61](#), [72](#)
- StopDacq
 - CMcsUsbDacqNet, [286](#)
- StopLoop
 - CMcsUsbDacqNet, [286](#)
- StopMCMovement
 - CMcsBus_MotorControlNet, [219](#)
- StopMovement
 - CRoboDeviceNet, [472](#)
- StopMovementI
 - CRoboStatorDeviceNet, [488](#)
- StopMovementXY
 - CRoboStatorDeviceNet, [488](#)
- StopMovementZ
 - CRoboStatorDeviceNet, [488](#)
- StopPlateClamp
 - CMultiwellDeviceNet, [381](#)
- StopPoll
 - CStimulusFunctionNet, [566](#)
 - CW2100_StimulatorFunctionNet, [610](#)
- StopSampling
 - CTEERFunctionNet, [593](#)
- StopTable
 - CRoboDacqNet, [459](#), [460](#)
- StorageCharge
 - Mcs::Usb, [67](#)
- StoreValveTable
 - CWarnerValveControllerDeviceNet, [651](#)
- SubtractFromAll
 - Mcs::Usb, [77](#)
- SubtractFromAllIOther
 - Mcs::Usb, [77](#)
- SubtractFromReferenceElectrodeOnly
 - Mcs::Usb, [77](#)
- SubtractionOff
 - Mcs::Usb, [77](#)
- SuperSpeed
 - Mcs::Usb, [69](#)
- Sw2to64
 - Mcs::Usb, [75](#)
- SwitchOnOff
 - CPositionIIDeviceNet, [409](#)
- SYNC_BIT0
 - CW2100_StimulatorFunctionNet, [611](#)
- SYNC_BIT1
 - CW2100_StimulatorFunctionNet, [611](#)
- SYNCOUT1
 - Mcs::Usb, [81](#)
- SYNCOUT2
 - Mcs::Usb, [81](#)
- SYNCOUT3
 - Mcs::Usb, [81](#)
- SYNCOUT4
 - Mcs::Usb, [81](#)
- SYNCOUT5
 - Mcs::Usb, [81](#)
- SYNCOUT6
 - Mcs::Usb, [81](#)
- SYNCOUT7
 - Mcs::Usb, [81](#)
- SYNCOUT8
 - Mcs::Usb, [81](#)
- syncoutdata
 - Mcs::Usb, [82](#)
- SyncStart
 - Mcs::Usb, [81](#)
- Table
 - Mcs::Usb, [90](#)
- Table_Wait
 - CRoboDacqNet, [460](#)
- TableDefBegin
 - CRoboDacqNet, [460](#)
- TableDefEnd
 - CRoboDacqNet, [460](#)
- TableEntryChangedEvent
 - CWarnerValveControllerDeviceNet, [654](#)
- TactSwitchGetState
 - CMcsBus_SensorNet, [229](#)
- TactSwitchSetDisplay
 - CMcsBus_SensorNet, [229](#)
- TBSI_127
 - Mcs::Usb, [50](#)
- TBSI_15
 - Mcs::Usb, [50](#)
- TBSI_31
 - Mcs::Usb, [50](#)
- TBSI_5
 - Mcs::Usb, [50](#)
- TBSI_63
 - Mcs::Usb, [50](#)
- TBSI_Dacq
 - Mcs::Usb, [75](#)
- TBSI_DACQDigitalSourceEnumNet
 - Mcs::Usb, [82](#)
- TBSI_Reserved
 - Mcs::Usb, [50](#)
- TbsiDacq
 - Mcs::Usb, [65](#)

- TbsiDacqHeadstage
 - Mcs::Usb, [64](#)
- TbsiDacqInterfaceboard
 - Mcs::Usb, [64](#)
- TC01
 - Mcs::Usb, [74](#)
- TC02
 - Mcs::Usb, [74](#)
- TCX
 - Mcs::Usb, [74](#)
- TcxDeviceTypeEnumNet
 - Mcs::Usb, [83](#)
- TcxSensorTypeEnumNet
 - Mcs::Usb, [83](#)
- TeerClampModeEnumNet
 - Mcs::Usb, [83](#)
- TEERFunctionNet
 - CTEERMachineDeviceNet, [594](#)
- TeerWaveformEnumNet
 - Mcs::Usb, [84](#)
- Tersens
 - Mcs::Usb, [74](#)
- Test_ADC_EPC10
 - Mcs::Usb, [66](#)
- Test_DAC_EPC10
 - Mcs::Usb, [66](#)
- ThrowCUsbExceptionNetOnError
 - CMcsUsbFunctionNet, [299](#)
 - CMcsUsbNet, [322](#)
- TiggerStatus
 - StgStatusNet, [695](#)
- TimeResolutionInNanoSeconds
 - W2100_StimulusParametersNet, [698](#)
- Timestamp
 - Mcs::Usb, [56](#)
- ToCpp
 - CFilterCoefficientsNet::s_FilterAttributesNet, [692](#)
- ToString
 - CFilterPropertyNet, [135](#)
 - CMcsUsbListEntryNet, [304](#)
 - HeadStageIDType, [678](#)
 - HeadStageIDTypeObject, [680](#)
- TouchTest
 - Mcs::Usb, [90](#)
- Triggerbox_AMS
 - Mcs::Usb, [74](#)
- Triggerbox_AMS3
 - Mcs::Usb, [74](#)
- Triggerbox_IMS
 - Mcs::Usb, [74](#)
- Triggerbox_R5
 - Mcs::Usb, [74](#)
- TriggerMask_Default
 - CRoboDacqNet, [460](#)
- TriggerOnly
 - Mcs::Usb, [81](#)
- TriggerSourceEnumNet
 - Mcs::Usb, [84](#)
- TriggerStatus
 - CMcsUsbDeviceStatePushFunctionNet, [288](#)
 - CMcsUsbDeviceStatePushNet, [289](#)
- TriggerStatus1
 - Mcs::Usb, [55](#)
- TriggerStatus2
 - Mcs::Usb, [55](#)
- TriggerStatus3
 - Mcs::Usb, [55](#)
- TriggerStatus4
 - Mcs::Usb, [55](#)
- TriggerValue_MoveAbs
 - CRoboDacqNet, [460](#)
- TriggerValue_StartQueue
 - CRoboDacqNet, [460](#)
- tsAuxIn1
 - Mcs::Usb, [85](#)
- tsAuxIn2
 - Mcs::Usb, [85](#)
- tsDACQCy1Dev1Runs
 - Mcs::Usb, [86](#)
- tsDACQCy1Dev2Runs
 - Mcs::Usb, [86](#)
- tsDACQCy2Dev1Runs
 - Mcs::Usb, [86](#)
- tsDACQCy2Dev2Runs
 - Mcs::Usb, [86](#)
- tsDigitalIn1
 - Mcs::Usb, [84](#)
- tsDigitalIn10
 - Mcs::Usb, [84](#)
- tsDigitalIn11
 - Mcs::Usb, [84](#)
- tsDigitalIn12
 - Mcs::Usb, [84](#)
- tsDigitalIn13
 - Mcs::Usb, [84](#)
- tsDigitalIn14
 - Mcs::Usb, [84](#)
- tsDigitalIn15
 - Mcs::Usb, [84](#)
- tsDigitalIn16
 - Mcs::Usb, [84](#)
- tsDigitalIn17
 - Mcs::Usb, [84](#)
- tsDigitalIn18
 - Mcs::Usb, [84](#)
- tsDigitalIn19
 - Mcs::Usb, [84](#)
- tsDigitalIn2
 - Mcs::Usb, [84](#)
- tsDigitalIn20
 - Mcs::Usb, [84](#)
- tsDigitalIn21
 - Mcs::Usb, [84](#)
- tsDigitalIn22
 - Mcs::Usb, [84](#)
- tsDigitalIn23
 - Mcs::Usb, [84](#)

Mcs::Usb, [84](#)
tsDigitalIn24
Mcs::Usb, [84](#)
tsDigitalIn25
Mcs::Usb, [85](#)
tsDigitalIn26
Mcs::Usb, [85](#)
tsDigitalIn27
Mcs::Usb, [85](#)
tsDigitalIn28
Mcs::Usb, [85](#)
tsDigitalIn29
Mcs::Usb, [85](#)
tsDigitalIn3
Mcs::Usb, [84](#)
tsDigitalIn30
Mcs::Usb, [85](#)
tsDigitalIn31
Mcs::Usb, [85](#)
tsDigitalIn32
Mcs::Usb, [85](#)
tsDigitalIn4
Mcs::Usb, [84](#)
tsDigitalIn5
Mcs::Usb, [84](#)
tsDigitalIn6
Mcs::Usb, [84](#)
tsDigitalIn7
Mcs::Usb, [84](#)
tsDigitalIn8
Mcs::Usb, [84](#)
tsDigitalIn9
Mcs::Usb, [84](#)
tsDigitalPuse0
Mcs::Usb, [85](#)
tsDigitalPuse1
Mcs::Usb, [85](#)
tsDigitalPuse10
Mcs::Usb, [86](#)
tsDigitalPuse11
Mcs::Usb, [86](#)
tsDigitalPuse12
Mcs::Usb, [86](#)
tsDigitalPuse13
Mcs::Usb, [86](#)
tsDigitalPuse14
Mcs::Usb, [86](#)
tsDigitalPuse15
Mcs::Usb, [86](#)
tsDigitalPuse16
Mcs::Usb, [86](#)
tsDigitalPuse17
Mcs::Usb, [86](#)
tsDigitalPuse18
Mcs::Usb, [86](#)
tsDigitalPuse19
Mcs::Usb, [86](#)
tsDigitalPuse2
Mcs::Usb, [85](#)
tsDigitalPuse20
Mcs::Usb, [86](#)
tsDigitalPuse21
Mcs::Usb, [86](#)
tsDigitalPuse22
Mcs::Usb, [86](#)
tsDigitalPuse23
Mcs::Usb, [86](#)
tsDigitalPuse24
Mcs::Usb, [86](#)
tsDigitalPuse25
Mcs::Usb, [86](#)
tsDigitalPuse26
Mcs::Usb, [86](#)
tsDigitalPuse27
Mcs::Usb, [86](#)
tsDigitalPuse28
Mcs::Usb, [86](#)
tsDigitalPuse29
Mcs::Usb, [86](#)
tsDigitalPuse3
Mcs::Usb, [85](#)
tsDigitalPuse30
Mcs::Usb, [86](#)
tsDigitalPuse31
Mcs::Usb, [86](#)
tsDigitalPuse4
Mcs::Usb, [85](#)
tsDigitalPuse5
Mcs::Usb, [85](#)
tsDigitalPuse6
Mcs::Usb, [85](#)
tsDigitalPuse7
Mcs::Usb, [85](#)
tsDigitalPuse8
Mcs::Usb, [86](#)
tsDigitalPuse9
Mcs::Usb, [86](#)
tsFeedback1
Mcs::Usb, [85](#)
tsFeedback10
Mcs::Usb, [85](#)
tsFeedback11
Mcs::Usb, [85](#)
tsFeedback12
Mcs::Usb, [85](#)
tsFeedback13
Mcs::Usb, [85](#)
tsFeedback14
Mcs::Usb, [85](#)
tsFeedback15
Mcs::Usb, [85](#)
tsFeedback16
Mcs::Usb, [85](#)
tsFeedback17
Mcs::Usb, [85](#)
tsFeedback18

- Mcs::Usb, [85](#)
- tsFeedback19
 - Mcs::Usb, [85](#)
- tsFeedback2
 - Mcs::Usb, [85](#)
- tsFeedback20
 - Mcs::Usb, [85](#)
- tsFeedback21
 - Mcs::Usb, [85](#)
- tsFeedback22
 - Mcs::Usb, [85](#)
- tsFeedback23
 - Mcs::Usb, [85](#)
- tsFeedback24
 - Mcs::Usb, [85](#)
- tsFeedback25
 - Mcs::Usb, [85](#)
- tsFeedback26
 - Mcs::Usb, [85](#)
- tsFeedback27
 - Mcs::Usb, [85](#)
- tsFeedback28
 - Mcs::Usb, [85](#)
- tsFeedback29
 - Mcs::Usb, [85](#)
- tsFeedback3
 - Mcs::Usb, [85](#)
- tsFeedback30
 - Mcs::Usb, [85](#)
- tsFeedback31
 - Mcs::Usb, [85](#)
- tsFeedback32
 - Mcs::Usb, [85](#)
- tsFeedback4
 - Mcs::Usb, [85](#)
- tsFeedback5
 - Mcs::Usb, [85](#)
- tsFeedback6
 - Mcs::Usb, [85](#)
- tsFeedback7
 - Mcs::Usb, [85](#)
- tsFeedback8
 - Mcs::Usb, [85](#)
- tsFeedback9
 - Mcs::Usb, [85](#)
- tsNone
 - Mcs::Usb, [84](#)
- tsSidebandBit8
 - Mcs::Usb, [86](#)
- tsTriggered
 - Mcs::Usb, [86](#)
- TxnGetSerialNumber
 - CMcsUsbNet, [322](#)
- TxnSetSerialNumber
 - CMcsUsbNet, [322](#)
- TxnTestMemoryReadAndCheck
 - CMcsUsbNet, [322](#)
- TxnTestMemoryWrite
 - CMcsUsbNet, [322](#)
- Type
 - HeadStageIDType, [679](#)
- TypeValue
 - HeadStageIDType, [679](#)
- Unknown
 - HeadStageIDType, [677](#)
 - Mcs::Usb, [50](#), [51](#), [63](#), [64](#), [83](#), [86](#)
- unknown
 - Mcs::Usb, [62](#)
- UnknownDest
 - Mcs::Usb, [54](#)
- UnknownSpeed
 - Mcs::Usb, [69](#)
- Unlock
 - Mcs::Usb, [73](#)
- UnlockPlateClamp
 - CMultiwellDeviceNet, [381](#)
- UnrolledAmplitude
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- UnrolledDuration
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- UnrolledSync
 - CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData, [696](#)
- Unsigned_16bit
 - Mcs::Usb, [58](#)
- Unsigned_24bit
 - Mcs::Usb, [58](#)
- Unsigned_32bit
 - Mcs::Usb, [58](#)
- UpdateChannelBlock
 - CCMOSMeaDeviceNet, [110](#)
- UpdateDisplay
 - CRoboDacqNet, [460](#)
- UpdateFirmware
 - CMcsUsbFactoryNet, [295–297](#)
- UpdateTransistorVoltages
 - CCMOSMea_FunctionNet, [107](#)
- UpdateTrigger
 - Mcs::Usb, [81](#)
- USB
 - FirmwareDestinationNames, [676](#)
 - Mcs::Usb, [51](#)
- USB_TARGET1
 - Mcs::Usb, [54](#)
- USB_TARGET2
 - Mcs::Usb, [54](#)
- USB_TARGET3
 - Mcs::Usb, [54](#)
- usbSetupPacket_t, [697](#)
 - bmRequestType, [697](#)
 - bRequest, [697](#)
 - wIndex, [697](#)
 - wLength, [697](#)
 - wValue, [697](#)

UsbTest
 Mcs::Usb, [75](#)
 UsbVendorIdEnumNet
 Mcs::Usb, [86](#)
 User_ADC_0
 Mcs::Usb, [66](#)
 User_ADC_1
 Mcs::Usb, [66](#)
 User_ADC_2
 Mcs::Usb, [66](#)
 User_ADC_3
 Mcs::Usb, [66](#)
 User_ADC_4
 Mcs::Usb, [66](#)
 User_DAC_0
 Mcs::Usb, [66](#)
 User_DAC_1
 Mcs::Usb, [66](#)
 User_DAC_2
 Mcs::Usb, [66](#)
 UserDefinedName
 HeadStageIDType, [679](#)
 UssingChamber
 Mcs::Usb, [65](#)
 UssingClampModeEnumNet
 Mcs::Usb, [86](#)
 UssingRail
 Mcs::Usb, [65](#)
 UssingUnitEnumNet
 Mcs::Usb, [87](#)

 Valid
 HeadStageIDType, [679](#)
 ValidKey
 CMcsUsbNet, [322](#)
 VendorIdEnumNet
 Mcs::Usb, [87](#)
 VendorInRequest
 CGenericDevelopDeviceNet, [162](#)
 VendorOutRequest
 CGenericDevelopDeviceNet, [162](#)
 VirtualDevice_ContinuousDacq
 CRoboDacqNet, [460](#)
 VirtualDevice_TableRun
 CRoboDacqNet, [460](#)
 Volt
 Mcs::Usb, [51](#), [87](#)
 Voltage
 BatteryState, [93](#)
 Voltage_3V3
 Mcs::Usb, [65](#)
 Voltage_5V0
 Mcs::Usb, [65](#)
 VoltageClamp
 Mcs::Usb, [87](#)
 VoltageRangeDisplayStringMilliVolt
 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet,
 [597](#)
 VoltageRangeInMicroVolt
 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet,
 [597](#)
 W2100_StimulusParametersNet, [698](#)
 VoltageResolutionInMicroVolt
 W2100_StimulusParametersNet, [698](#)
 VoltageString
 BatteryState, [93](#)
 VOPSTimerSetResetTimes
 CCMOSMea_FunctionNet, [107](#)

 W16IsW14
 HeadStageIDType, [679](#)
 W2100
 Mcs::Usb, [75](#)
 W2100_Accel_Gyro_Select_EnumNet
 Mcs::Usb, [87](#)
 W2100_FunctionNet
 CMeaDeviceNet, [351](#)
 W2100_StimulusParametersNet, [697](#)
 CurrentRangeInNanoAmp, [698](#)
 CurrentResolutionInNanoAmp, [698](#)
 DACResolution, [698](#)
 TimeResolutionInNanoSeconds, [698](#)
 VoltageRangeInMicroVolt, [698](#)
 VoltageResolutionInMicroVolt, [698](#)
 W2100DacqGroupChannelEnumNet
 Mcs::Usb, [87](#)
 W2100DigitalSourceEnumNet
 Mcs::Usb, [89](#)
 W2100IFB2
 Mcs::Usb, [65](#)
 W2100Interfaceboard
 Mcs::Usb, [64](#)
 W2100WirelessReceiver
 Mcs::Usb, [64](#), [65](#)
 W2100WirelessReceiverAnalog
 Mcs::Usb, [64](#), [65](#)
 WaitForAllChambers
 CWarnerUssingFunctionNet, [629](#)
 WaitForChamber
 CWarnerUssingFunctionNet, [629](#)
 WaitTimer
 CRoboDeviceNet, [472](#)
 Warner
 Mcs::Usb, [83](#)
 Warner_TEER_Machine
 Mcs::Usb, [76](#)
 Warner_Ussing
 Mcs::Usb, [76](#)
 WARNER_USSING_DEVICE
 Mcs::Usb, [59](#)
 Warner_Valve_Control
 Mcs::Usb, [76](#)
 WARNER_VALVE_CONTROL_DEVICE
 Mcs::Usb, [59](#)
 WarnerUssingFunction
 CWarnerUssingDeviceNet, [613](#)
 WClassicFunctionNet
 CMeaDeviceNet, [351](#)

Whole_Cell_Patch
Mcs::Usb, [75](#)

WholeCellPatch
Mcs::Usb, [78](#)

WholeCellPatchHeadstage
Mcs::Usb, [65](#)

wIndex
usbSetupPacket_t, [697](#)

WirelessHeadStageAccDataRE1HS1
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE1HS2
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE1HS3
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE1HS4
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE2HS1
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE2HS2
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE2HS3
Mcs::Usb, [88](#)

WirelessHeadStageAccDataRE2HS4
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE1HS1
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE1HS2
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE1HS3
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE1HS4
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE2HS1
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE2HS2
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE2HS3
Mcs::Usb, [88](#)

WirelessHeadStageAnalogRE2HS4
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE1HS1
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE1HS2
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE1HS3
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE1HS4
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE2HS1
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE2HS2
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE2HS3
Mcs::Usb, [88](#)

WirelessHeadStageGyroDataRE2HS4
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE1HS1
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE1HS2
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE1HS3
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE1HS4
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE2HS1
Mcs::Usb, [88](#)

WirelessHeadStageOptoStimCurrentRE2HS2
Mcs::Usb, [89](#)

WirelessHeadStageOptoStimCurrentRE2HS3
Mcs::Usb, [89](#)

WirelessHeadStageOptoStimCurrentRE2HS4
Mcs::Usb, [89](#)

WirelessHeadStageReservedARE1HS1
Mcs::Usb, [88](#)

WirelessHeadStageReservedARE1HS2
Mcs::Usb, [88](#)

WirelessHeadStageReservedARE1HS3
Mcs::Usb, [88](#)

WirelessHeadStageReservedARE1HS4
Mcs::Usb, [88](#)

WirelessHeadStageReservedARE2HS1
Mcs::Usb, [89](#)

WirelessHeadStageReservedARE2HS2
Mcs::Usb, [89](#)

WirelessHeadStageReservedARE2HS3
Mcs::Usb, [89](#)

WirelessHeadStageReservedARE2HS4
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE1HS1
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE1HS2
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE1HS3
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE1HS4
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE2HS1
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE2HS2
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE2HS3
Mcs::Usb, [89](#)

WirelessHeadStageReservedBRE2HS4
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE1HS1
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE1HS2
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE1HS3
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE1HS4
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE2HS1
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE2HS2
Mcs::Usb, [89](#)

WirelessHeadStageReservedCRE2HS3
 Mcs::Usb, [89](#)
 WirelessHeadStageReservedCRE2HS4
 Mcs::Usb, [89](#)
 WirelessHeadStageStatusRE1HS1
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE1HS2
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE1HS3
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE1HS4
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE2HS1
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE2HS2
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE2HS3
 Mcs::Usb, [88](#)
 WirelessHeadStageStatusRE2HS4
 Mcs::Usb, [88](#)
 WirelessTestAdapter
 Mcs::Usb, [50](#)
 wLength
 usbSetupPacket_t, [697](#)
 Work
 Mcs::Usb, [90](#)
 WPA16
 Mcs::Usb, [75](#)
 WPA32
 Mcs::Usb, [75](#)
 WPA4
 Mcs::Usb, [75](#)
 WPA8
 Mcs::Usb, [75](#)
 WPAError_ScanningIsPending
 CMcsUsbNet, [328](#)
 Write
 CExternDTesterDeviceNet, [128](#)
 Write2
 CExternDTesterDeviceNet, [129](#)
 WriteEepromRegisterPreconfig
 CMcsUsbNet, [323](#)
 WritePipe
 CGenericDevelopDeviceNet, [162](#)
 WriteRegister
 CMcsUsbNet, [323](#)
 WriteRegister32
 CMcsUsbNet, [323](#)
 WriteRegisterArray
 CMcsUsbNet, [323](#)
 WriteRegisterTimeSlot
 CMcsUsbNet, [324](#)
 WriteRegisterValue
 CMcsUsbNet, [324](#)
 WriteUARTData
 CLIH3DeviceNet, [197](#)
 wValue
 usbSetupPacket_t, [697](#)
 WvcDisplayModeEnumNet
 Mcs::Usb, [90](#)
 WvcValveModeEnumNet
 Mcs::Usb, [90](#)
 Zero
 Mcs::Usb, [60](#), [70](#), [79](#), [82](#), [89](#)