

McsUsbNet.dll

Version 5.1.7

Multi Channel Systems MCS GmbH
Aspenhaustrasse 21
72770 Reutlingen
Germany
Fon +49-71 21-90 92 5 - 0
Fax +49-71 21-90 92 5 -11
info@multichannelsystems.com
www.multichannelsystems.com

Generated by Doxygen 1.8.18

1	McsUsbNet.dll for MCS USB devices	1
	1.1 Introduction	1
	1.2 System requirements	2
	1.3 Connecting to an MCS device	2
2	Device Classes	2
	2.1 The MCS FluidControl Device	3
	2.1.1 Introduction	3
	2.1.2 Access to the FluidControl device	3
	2.2 MCS-USB-Sw2to64 device	3
3	Function Classes	4
4	Data ACQuisition (DACQ) Devices	5
5	The MCS Robo Device	6
	5.1 Introduction	6
6	STG200x & STG400x STimulus Generator	6
	6.1 Introduction	6
	6.2 Download mode	6
	6.2.1 Memory Layout and Trigger Setup	7
	6.3 Streaming mode	8
	6.3.1 Memory Layout and Trigger Setup	9
7	Namespace Index	11
	7.1 Namespace List	11
8	Hierarchical Index	11
	8.1 Class Hierarchy	11
9	Class Index	16
	9.1 Class List	16
1	0 Namespace Documentation	22
	10.1 Mcs Namespace Reference	22
	10.2 Mcs::Usb Namespace Reference	22
	10.2.1 Enumeration Type Documentation	27
	10.2.2 Function Documentation	27
1	1 Class Documentation	29
	11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference	29
	11.1.1 Member Data Documentation	29
	11.2 BatteryState Class Reference	29
	11.2.1 Property Documentation	30
	11.3 BesselFilterHighPassNet Class Reference	30

11.3.1 Constructor & Destructor Documentation	30
11.4 BesselFilterLowPassNet Class Reference	31
11.4.1 Constructor & Destructor Documentation	31
11.5 ButterworthFilterHighPassNet Class Reference	31
11.5.1 Constructor & Destructor Documentation	32
11.6 ButterworthFilterLowPassNet Class Reference	32
11.6.1 Constructor & Destructor Documentation	32
11.7 CChannelTestDeviceNet Class Reference	33
11.7.1 Constructor & Destructor Documentation	33
11.7.2 Member Function Documentation	33
11.8 CCMOSMea_FunctionNet Class Reference	34
11.8.1 Constructor & Destructor Documentation	36
11.8.2 Member Function Documentation	36
11.9 CCMOSMeaDeviceNet Class Reference	45
11.9.1 Constructor & Destructor Documentation	46
11.9.2 Member Function Documentation	46
11.9.3 Property Documentation	47
11.10 CCreateFilterNet Class Reference	48
11.10.1 Constructor & Destructor Documentation	48
11.10.2 Member Function Documentation	49
11.10.3 Property Documentation	49
11.11 CDacCalibrationFunctionNet Class Reference	50
11.11.1 Detailed Description	50
11.11.2 Constructor & Destructor Documentation	50
11.11.3 Member Function Documentation	51
11.12 CDacqGroupChannelGenericSelectionNet Class Reference	52
11.12.1 Constructor & Destructor Documentation	52
11.13 CDacqGroupChannelSelectionNet Class Reference	52
11.13.1 Constructor & Destructor Documentation	53
$11.14 CDacqGroupChannelSelectionTemplateNet < DacqGroupChannelEnumTemplateNet, Dacq \leftarrow$	
GroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference	53
11.14.1 Constructor & Destructor Documentation	54
11.14.2 Member Function Documentation	54
11.15 CDeviceGroupChannelInfoGenericNet Class Reference	56
11.15.1 Constructor & Destructor Documentation	56
11.16 CDeviceGroupChannelInfoNet Class Reference	56
11.16.1 Constructor & Destructor Documentation	57
11.17 CDeviceGroupChannelInfoSCUNet Class Reference	57
11.17.1 Constructor & Destructor Documentation	57
11.18 CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet > Class Template Reference	57
11.18.1 Constructor & Destructor Documentation	58
11.18.2 Member Data Documentation	58

11.19 CDeviceGroupChannelInfoW2100Net Class Reference	58
11.19.1 Constructor & Destructor Documentation	59
11.20 CDigOutStimulatorFunctionNet Class Reference	59
11.20.1 Detailed Description	60
11.20.2 Constructor & Destructor Documentation	60
11.20.3 Member Function Documentation	60
11.21 CEncapsulatorDeviceNet Class Reference	63
11.21.1 Detailed Description	63
11.21.2 Constructor & Destructor Documentation	64
11.21.3 Member Function Documentation	64
11.22 CExternDTesterDeviceNet Class Reference	64
11.22.1 Detailed Description	64
11.22.2 Constructor & Destructor Documentation	65
11.22.3 Member Function Documentation	65
11.23 CFilterCoefficientsNet Class Reference	66
11.23.1 Constructor & Destructor Documentation	66
11.23.2 Member Function Documentation	67
11.23.3 Property Documentation	68
11.24 CFilterConfigurationNet Class Reference	68
11.24.1 Constructor & Destructor Documentation	68
11.24.2 Member Function Documentation	69
11.25 CFilterConfigurationRegisterNet Class Reference	70
11.25.1 Constructor & Destructor Documentation	70
11.25.2 Member Function Documentation	70
11.26 CFilterPropertyNet Class Reference	72
11.26.1 Constructor & Destructor Documentation	72
11.26.2 Member Function Documentation	72
11.26.3 Property Documentation	72
11.27 CFluidControlDeviceNet Class Reference	73
11.27.1 Detailed Description	74
11.27.2 Constructor & Destructor Documentation	75
11.27.3 Member Function Documentation	75
11.27.4 Property Documentation	80
11.28 CFYIDeviceNet Class Reference	80
11.28.1 Detailed Description	81
11.28.2 Constructor & Destructor Documentation	81
11.28.3 Property Documentation	81
11.29 CGenericDevelopDeviceNet Class Reference	82
11.29.1 Detailed Description	88
11.29.2 Constructor & Destructor Documentation	89
11.29.3 Member Function Documentation	89
11.30 CGilson DaviceNet Class Reference	100

	11.30.1 Detailed Description				 	 	 	 	. 101
	11.30.2 Constructor & Destructor Documentation				 	 	 	 	. 101
	11.30.3 Member Function Documentation				 	 	 	 	. 101
	11.30.4 Member Data Documentation				 	 	 	 	. 102
11.3	CGrapheneFunctionNet Class Reference				 	 	 	 	. 102
	11.31.1 Detailed Description				 	 	 	 	. 103
	11.31.2 Constructor & Destructor Documentation				 	 	 	 	. 103
	11.31.3 Member Function Documentation				 	 	 	 	. 103
11.32	2 CHiClampDeviceNet Class Reference				 	 	 	 	. 109
	11.32.1 Detailed Description				 	 	 	 	. 110
	11.32.2 Constructor & Destructor Documentation				 	 	 	 	. 110
	11.32.3 Property Documentation				 	 	 	 	. 110
11.33	3 CHLADacqNet Class Reference				 	 	 	 	. 110
	11.33.1 Constructor & Destructor Documentation				 	 	 	 	. 111
11.34	CHLADeviceNet Class Reference				 	 	 	 	. 111
	11.34.1 Detailed Description				 	 	 	 	. 111
	11.34.2 Constructor & Destructor Documentation				 	 	 	 	. 111
	11.34.3 Property Documentation				 	 	 	 	. 112
11.35	5 CMcsUsbDacqNet::CHWInfo Class Reference				 	 	 	 	. 112
	11.35.1 Detailed Description				 	 	 	 	. 112
	11.35.2 Constructor & Destructor Documentation				 	 	 	 	. 113
	11.35.3 Member Function Documentation				 	 	 	 	. 113
11.36	6 CIntanMea_FunctionNet Class Reference				 	 	 	 	. 114
	11.36.1 Constructor & Destructor Documentation				 	 	 	 	. 115
	11.36.2 Member Function Documentation				 	 	 	 	. 115
11.37	7 CInterfaceboard2FunctionNet Class Reference .				 	 	 	 	. 116
	11.37.1 Detailed Description				 	 	 	 	. 117
	11.37.2 Constructor & Destructor Documentation				 	 	 	 	. 117
	11.37.3 Member Function Documentation				 	 	 	 	. 117
11.38	3 CInterfaceboardFunctionNet Class Reference				 	 	 	 	. 118
	11.38.1 Detailed Description				 	 	 	 	. 119
	11.38.2 Constructor & Destructor Documentation				 	 	 	 	. 119
	11.38.3 Member Function Documentation				 	 	 	 	. 119
11.39	OCLIH3DeviceNet Class Reference				 	 	 	 	. 120
	11.39.1 Detailed Description				 	 	 	 	. 122
	11.39.2 Constructor & Destructor Documentation				 	 	 	 	. 122
	11.39.3 Member Function Documentation				 	 	 	 	. 122
	11.39.4 Property Documentation				 	 	 	 	. 129
11.40	CMcsBus_AxisParametersNet Class Reference .				 	 	 	 	. 129
	11.40.1 Constructor & Destructor Documentation				 	 	 	 	. 129
	11.40.2 Member Function Documentation				 	 	 	 	. 130
11 4	CMcsBus ExtensionNet Class Reference				 	 	 	 	. 130

11.41.1 Constructor & Destructor Documentation	31
11.41.2 Member Function Documentation	31
11.42 CMcsBus_FYIExtensionNet Class Reference	31
11.42.1 Constructor & Destructor Documentation	32
11.42.2 Member Function Documentation	32
11.43 CMcsBus_MotorControlNet Class Reference	3
11.43.1 Constructor & Destructor Documentation	36
11.43.2 Member Function Documentation	36
11.44 CMcsBus_SensorNet Class Reference	50
11.44.1 Constructor & Destructor Documentation	52
11.44.2 Member Function Documentation	52
11.45 CMcsBus_TempSensorNet Class Reference	31
11.45.1 Constructor & Destructor Documentation	31
11.45.2 Member Function Documentation	32
11.46 CMcsBus_VoltageModeNet Class Reference	3
11.46.1 Constructor & Destructor Documentation	34
11.46.2 Member Function Documentation	34
11.47 CMcsBusNet Class Reference	67
11.47.1 Constructor & Destructor Documentation	8
11.47.2 Member Function Documentation	8
11.48 CMcsUsbDacqNet Class Reference	1
11.48.1 Detailed Description	7
11.48.2 Constructor & Destructor Documentation	7
11.48.3 Member Function Documentation	7
11.48.4 Member Data Documentation	8
11.48.5 Property Documentation	9
11.48.6 Event Documentation	9
11.49 CMcsUsbDeviceStatePushFunctionNet Class Reference	9
11.49.1 Constructor & Destructor Documentation	20
11.49.2 Member Function Documentation	20
11.49.3 Event Documentation	20
11.50 CMcsUsbDeviceStatePushNet Class Reference	20
11.50.1 Constructor & Destructor Documentation	21
11.50.2 Member Function Documentation	21
11.50.3 Event Documentation	21
11.51 CMcsUsbFactoryNet Class Reference	21
11.51.1 Constructor & Destructor Documentation	23
11.51.2 Member Function Documentation	23
11.51.3 Member Data Documentation	29
11.52 CMcsUsbFunctionNet Class Reference	30
11.52.1 Constructor & Destructor Documentation	30
11.52.2 Member Function Documentation	31

11.52.3 Member Data Documentation
11.53 CMcsUsbFunctionPointerContainer Class Reference
11.54 CMcsUsbListEntryNet Class Reference
11.54.1 Detailed Description
11.54.2 Constructor & Destructor Documentation
11.54.3 Member Function Documentation
11.54.4 Property Documentation
11.55 CMcsUsbListNet Class Reference
11.55.1 Detailed Description
11.55.2 Constructor & Destructor Documentation
11.55.3 Member Function Documentation
11.55.4 Property Documentation
11.55.5 Event Documentation
11.56 CMcsUsbNet Class Reference
11.56.1 Detailed Description
11.56.2 Constructor & Destructor Documentation
11.56.3 Member Function Documentation
11.56.4 Member Data Documentation
11.56.5 Property Documentation
11.57 CMcsUsbPointerContainer Class Reference
11.58 CMEA2100x256FunctionNet Class Reference
11.58.1 Detailed Description
11.58.2 Constructor & Destructor Documentation
11.58.3 Member Function Documentation
11.59 CMeaAudioFunctionNet Class Reference
11.59.1 Constructor & Destructor Documentation
11.59.2 Member Function Documentation
11.60 CMeaCleanDeviceNet Class Reference
11.60.1 Detailed Description
11.60.2 Constructor & Destructor Documentation
11.60.3 Member Function Documentation
11.61 CMeaCoatDeviceNet Class Reference
11.61.1 Detailed Description
11.61.2 Constructor & Destructor Documentation
11.61.3 Member Function Documentation
11.62 CMeaDeviceNet Class Reference
11.62.1 Detailed Description
11.62.2 Constructor & Destructor Documentation
11.62.3 Member Function Documentation
11.62.4 Property Documentation
11.63 CMeaDigitalDataFunctionNet Class Reference
11.63.1 Constructor & Destructor Documentation

11.63.2 Member Function Documentation
11.64 CMeaFeedbackFunctionNet Class Reference
11.64.1 Constructor & Destructor Documentation
11.64.2 Member Function Documentation
11.65 CMealmpedanceDeviceNet Class Reference
11.65.1 Constructor & Destructor Documentation
11.65.2 Member Function Documentation
11.66 CMeasureTableDeviceNet Class Reference
11.66.1 Detailed Description
11.66.2 Constructor & Destructor Documentation
11.66.3 Property Documentation
11.67 CMeaSwitchDeviceNet Class Reference
11.67.1 Detailed Description
11.67.2 Constructor & Destructor Documentation
11.67.3 Member Function Documentation
11.68 CMeaUSBDeviceNet Class Reference
11.68.1 Detailed Description
11.68.2 Constructor & Destructor Documentation
11.69 CMeFunctionNet Class Reference
11.69.1 Detailed Description
11.69.2 Constructor & Destructor Documentation
11.69.3 Member Function Documentation
11.70 CMultiBatteryChargerDeviceNet Class Reference
11.70.1 Detailed Description
11.70.2 Constructor & Destructor Documentation
11.70.3 Member Function Documentation
11.71 CMultiwellCallbackFunctionNet Class Reference
11.71.1 Detailed Description
11.71.2 Constructor & Destructor Documentation
11.71.3 Member Function Documentation
11.71.4 Event Documentation
11.72 CMultiwellDeviceNet Class Reference
11.72.1 Detailed Description
11.72.2 Constructor & Destructor Documentation
11.72.3 Member Function Documentation
11.73 CMultiwellOptoStimFunctionNet Class Reference
11.73.1 Detailed Description
11.73.2 Constructor & Destructor Documentation
11.73.3 Member Function Documentation
11.74 CNF_GenDeviceNet Class Reference
11.74.1 Constructor & Destructor Documentation
11.74.2 Member Function Documentation

11.75 COctoPotDeviceNet Class Reference	318
11.75.1 Constructor & Destructor Documentation	318
11.75.2 Member Function Documentation	319
11.76 COkuvisionStimulatorDeviceNet Class Reference	322
11.76.1 Constructor & Destructor Documentation	323
11.76.2 Member Function Documentation	323
11.77 CPatchServerDeviceNet Class Reference	325
11.77.1 Detailed Description	326
11.77.2 Constructor & Destructor Documentation	326
11.77.3 Property Documentation	326
11.78 CPathIdentDeviceNet Class Reference	326
11.78.1 Constructor & Destructor Documentation	327
11.78.2 Member Function Documentation	327
11.79 CPedoterDeviceNet Class Reference	327
11.79.1 Detailed Description	328
11.79.2 Constructor & Destructor Documentation	328
11.79.3 Member Function Documentation	328
11.80 CPeristalticPumpDeviceNet Class Reference	329
11.80.1 Detailed Description	329
11.80.2 Constructor & Destructor Documentation	329
11.80.3 Property Documentation	330
11.81 CPgaDeviceNet Class Reference	330
11.81.1 Constructor & Destructor Documentation	330
11.81.2 Member Function Documentation	331
11.82 CPositionIIDeviceNet Class Reference	332
11.82.1 Detailed Description	334
11.82.2 Constructor & Destructor Documentation	334
11.82.3 Member Function Documentation	334
11.82.4 Property Documentation	341
11.83 CPositionImpDeviceNet Class Reference	341
11.83.1 Detailed Description	341
11.83.2 Constructor & Destructor Documentation	342
11.83.3 Member Function Documentation	342
11.84 CPPCDeviceNet Class Reference	344
11.84.1 Constructor & Destructor Documentation	344
11.84.2 Property Documentation	345
11.85 CPPCFunctionNet Class Reference	345
11.85.1 Detailed Description	346
11.85.2 Constructor & Destructor Documentation	346
11.85.3 Member Function Documentation	347
11.86 CPPS_DeviceNet Class Reference	354
11.86.1 Constructor & Destructor Documentation	354

11.86.2 Property Documentation
11.87 CPPS_FunctionNet Class Reference
11.87.1 Constructor & Destructor Documentation
11.87.2 Member Function Documentation
11.88 CPPSDeviceNet Class Reference
11.88.1 Detailed Description
11.88.2 Constructor & Destructor Documentation
11.89 CProgramPressureCurveNet Class Reference
11.89.1 Detailed Description
11.89.2 Constructor & Destructor Documentation
11.89.3 Member Function Documentation
11.90 CPulseGeneratorFunctionNet Class Reference
11.90.1 Detailed Description
11.90.2 Constructor & Destructor Documentation
11.90.3 Member Function Documentation
11.91 CRadioControledDevicesNet Class Reference
11.91.1 Constructor & Destructor Documentation
11.91.2 Member Function Documentation
11.92 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference
11.92.1 Constructor & Destructor Documentation
11.92.2 Member Function Documentation
11.92.3 Member Data Documentation
11.93 CRetinaLedDeviceNet Class Reference
11.93.1 Constructor & Destructor Documentation
11.93.2 Member Function Documentation
11.94 CRFFunctionNet Class Reference
11.94.1 Detailed Description
11.94.2 Constructor & Destructor Documentation
11.94.3 Member Function Documentation
11.95 CRobo_FYIProgram_FunctionNet Class Reference
11.95.1 Constructor & Destructor Documentation
11.95.2 Member Function Documentation
11.96 CRobo_FYITemp_FunctionNet Class Reference
11.96.1 Constructor & Destructor Documentation
11.96.2 Member Function Documentation
11.97 CRoboDacqNet Class Reference
11.97.1 Constructor & Destructor Documentation
11.97.2 Member Function Documentation
11.98 CRoboDeviceNet Class Reference
11.98.1 Detailed Description
11.98.2 Constructor & Destructor Documentation
11.98.3 Member Function Documentation 39

11.98.4 Member Data Documentation
11.98.5 Property Documentation
11.98.6 Event Documentation
11.99 CRoboFluidDeviceNet Class Reference
11.99.1 Constructor & Destructor Documentation
11.99.2 Member Function Documentation
11.99.3 Member Data Documentation
11.99.4 Property Documentation
11.100 CRoboInjectDeviceNet Class Reference
11.100.1 Detailed Description
11.100.2 Constructor & Destructor Documentation
11.101 CRoboocyte2DeviceNet Class Reference
11.101.1 Detailed Description
11.101.2 Constructor & Destructor Documentation
11.101.3 Member Function Documentation
11.102 CRoboStatorDeviceNet Class Reference
11.102.1 Constructor & Destructor Documentation
11.102.2 Member Function Documentation
11.102.3 Property Documentation
11.103 CSafeISDeviceNet Class Reference
11.103.1 Detailed Description
11.103.2 Constructor & Destructor Documentation
11.103.3 Member Function Documentation
11.103.4 Property Documentation
11.104 CSCUDacqGroupChannelSelectionNet Class Reference
11.104.1 Constructor & Destructor Documentation
11.105 CSCUFunctionNet Class Reference
11.105.1 Detailed Description
11.105.2 Constructor & Destructor Documentation
11.105.3 Member Function Documentation
11.105.4 Event Documentation
11.106 CSerialPortNet Class Reference
11.106.1 Constructor & Destructor Documentation
11.106.2 Member Function Documentation
11.107 CStg200xBasicNet Class Reference
11.107.1 Detailed Description
11.107.2 Constructor & Destructor Documentation
11.107.3 Member Function Documentation
11.108 CStg200xDownloadBasicNet Class Reference
11.108.1 Detailed Description
11.108.2 Member Function Documentation
11.108.3 Property Documentation 47

11.109 CStg200xDownloadNet Class Reference
11.109.1 Detailed Description
11.109.2 Constructor & Destructor Documentation
11.109.3 Member Function Documentation
11.109.4 Event Documentation
11.110 CStimulusFunctionNet Class Reference
11.110.1 Constructor & Destructor Documentation
11.110.2 Member Function Documentation
11.110.3 Event Documentation
11.111 CSw2to64DeviceNet Class Reference
11.111.1 Detailed Description
11.111.2 Constructor & Destructor Documentation
11.111.3 Member Function Documentation
11.112 CTcxDeviceNet Class Reference
11.112.1 Detailed Description
11.112.2 Constructor & Destructor Documentation
11.112.3 Member Function Documentation
11.113 CTEERFunctionNet Class Reference
11.113.1 Detailed Description
11.113.2 Constructor & Destructor Documentation
11.113.3 Member Function Documentation
11.114 CTEERMachineDeviceNet Class Reference
11.114.1 Constructor & Destructor Documentation
11.114.2 Property Documentation
11.115 CUsbExceptionNet Class Reference
11.115.1 Detailed Description
11.115.2 Constructor & Destructor Documentation
11.115.3 Property Documentation
11.116 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference
11.116.1 Constructor & Destructor Documentation
11.116.2 Member Data Documentation
11.117 CW2100_FunctionNet Class Reference
11.117.1 Constructor & Destructor Documentation
11.117.2 Member Function Documentation
11.117.3 Property Documentation
11.118 CW2100_StimulatorFunctionNet Class Reference
11.118.1 Constructor & Destructor Documentation
11.118.2 Member Function Documentation
11.118.3 Member Data Documentation
11.118.4 Event Documentation
11.119 CW2100DacqGroupChannelSelectionNet Class Reference
11.119.1 Constructor & Destructor Documentation

11.120 CWarnerUssingDeviceNet Class Reference	40
11.120.1 Detailed Description	40
11.120.2 Constructor & Destructor Documentation	40
11.120.3 Property Documentation	41
11.121 CWarnerUssingFunctionNet Class Reference	41
11.121.1 Detailed Description	43
11.121.2 Constructor & Destructor Documentation	43
11.121.3 Member Function Documentation	43
11.122 CWarnerValveControllerDeviceNet Class Reference	56
11.122.1 Detailed Description	60
11.122.2 Constructor & Destructor Documentation	60
11.122.3 Member Function Documentation	61
11.122.4 Event Documentation	77
11.123 CWarnerValveControllerDeviceTesterFunctionNet Class Reference	80
11.123.1 Detailed Description	81
11.123.2 Constructor & Destructor Documentation	81
11.123.3 Member Function Documentation	81
11.124 CWClassicFunctionNet Class Reference	83
11.124.1 Constructor & Destructor Documentation	84
11.124.2 Member Function Documentation	84
11.125 CWirelessBaseFunctionNet Class Reference	88
11.125.1 Constructor & Destructor Documentation	88
11.125.2 Member Function Documentation	88
11.126 DeviceIdNet Struct Reference	88
11.126.1 Detailed Description	89
11.126.2 Constructor & Destructor Documentation	89
11.126.3 Member Function Documentation	89
11.126.4 Member Data Documentation	90
11.127 DigitalSource< digitalsourceenum > Class Template Reference	90
11.127.1 Constructor & Destructor Documentation	90
11.127.2 Member Function Documentation	90
11.127.3 Property Documentation	91
11.128 DriverVersionNet Class Reference	91
11.128.1 Detailed Description	92
11.128.2 Constructor & Destructor Documentation	92
11.128.3 Member Function Documentation	92
11.129 FirmwareDestinationNames Class Reference	96
11.129.1 Member Data Documentation	97
11.130 HeadStageIDType Class Reference	00
11.130.1 Member Enumeration Documentation	00
11.130.2 Constructor & Destructor Documentation	01
11 130 3 Member Function Documentation	Λ1

622

11.130.4 Property Documentation
11.131 HeadstageIDTypeObject Class Reference
11.131.1 Constructor & Destructor Documentation
11.131.2 Member Function Documentation
11.131.3 Member Data Documentation
11.131.4 Property Documentation
11.132 HeadStageIDTypeState Class Reference
11.132.1 Property Documentation
11.133 mkfilterNet Class Reference
11.133.1 Member Function Documentation
11.134 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference 608
11.134.1 Member Function Documentation
11.135 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference 615
11.135.1 Member Function Documentation
11.136 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference 615
11.136.1 Constructor & Destructor Documentation
11.136.2 Member Function Documentation
11.136.3 Member Data Documentation
11.137 CMeaAudioFunctionNet::s_setaudionet Struct Reference 616
11.137.1 Member Data Documentation
11.138 CStimulusFunctionNet::SidebandData Class Reference 617
11.138.1 Constructor & Destructor Documentation
11.138.2 Property Documentation
11.139 StgStatusNet Class Reference
11.139.1 Member Function Documentation
11.139.2 Member Data Documentation
11.140 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference 618
11.140.1 Constructor & Destructor Documentation
11.140.2 Property Documentation
11.141 usbSetupPacket_t Class Reference
11.141.1 Member Data Documentation
11.142 W2100_StimulusParametersNet Struct Reference
11.142.1 Member Data Documentation

1 McsUsbNet.dll for MCS USB devices

1.1 Introduction

Index

This DLL provides the .NET interface to MCS devices

The most important options are accessing our stimulator and data acquisition devices:

- STG200x & STG400x STimulus Generator
- · Data ACQuisition (DACQ) Devices

See here for a list of our other devices: Device Classes.

And here for a list of function classes addressing groups of features that might be shared between different devices: Function Classes.

1.2 System requirements

The DLL can be used with any .NET compatible language.

The DLL needs the .NET Framework 4.7.2.

It requires the Microsoft Visual C++ Redistributable for Visual Studio 2019 to be installed.

It also requires the **USB driver** to be installed.

The simplest way to achieve this is to install the latest **Multi Channel Experimenter** setup (will install 64bit redistributable).

All examples assume that the Mcs.Usb namespace is loaded: using namespace Mcs.Usb;

Include the file McsUsbNet.dll into the references of your project.

1.3 Connecting to an MCS device

A connection to a DAQ device is established by Mcs.Usb.CMcsUsbNet.Connect. When this function is called without argument, the first DAQ device found on the USB bus is used:

```
CMcsUsbNet device = new CMcsUsbNet();
device.Connect();
```

When more than one DAQ device of the specific type is connected, you can use the Mcs.Usb.CMcsUsbListNet class to get a list of available devices:

```
CMcsUsbListNet usblist = new CMcsUsbListNet(DeviceEnumNet.MCS_DEVICE_USB);
var entry = usblist.GetUsbListEntry((uint)0);
CMcsUsbNet device = new CMcsUsbNet();
device.Connect(entry);
```

After you are finished with the device, you can disconnect the device object from the device by: device.Disonnect();

2 Device Classes

- For FluidControl device see MCS FluidControl
- For SW2TO64 device see MCS-USB-Sw2to64
- For TCx device see Mcs.Usb.CTcxDeviceNet

2.1 The MCS FluidControl Device

2.1.1 Introduction

The FluidControl Device can control up to 24 valves. The nominal voltage is 24V.

8 TTL level digital output ports are available and 8 TTL inputs can be read in.

The device has 8 ADC inputs with a rage from 0V to 3.3V.

2.1.2 Access to the FluidControl device

For connecting to a FluidControl device see Connecting to an MCS device.*

CFluidControlDevice* m_dacq;

m_fluidcontrol = new CFluidControlDevice;

status = m fluidcontrol->Connect();

The valves are controlled with the CFluidControlDevice::SetValve call. The argument given is a bit pattern of all valves which should be open.

The digital outputs can be controlled with the CFluidControlDevice::SetDigout call. Again, a bit pattern of all digital output pins which should be set to a logic high level is given as an argument.

The current state of the valves and the digital outputs can be read back with the CFluidControlDevice::GetValve and CFluidControlDevice::GetDigout

The command to read an ADC-Channel is CFluidControlDevice::GetAdc. Here the channelnummer which should be read in is given as an argument and the return value is the current Adc level.

The state of the digital inputs is read with the CFluidControlDevice::GetDigin call. Here the return value is the bit pattern of the digital inputs.

The connection to the device is closed with the CFluidControlDevice::Disconnect call.

2.2 MCS-USB-Sw2to64 device

The class Mcs.Usb.CSw2to64DeviceNet controls the setting of the switches in the MCS-USB-Sw2to64 device.

First construct an object of the class:

```
CSw2to64DeviceNet device = new CSw2to64DeviceNet();
```

For connecting to an MCS-USB-Sw2to64 device see Connecting to an MCS device.

To get the number of channels the device handles:

```
int number = device.GetNumber()
```

Set all channel switches at once:

```
byte[] pattern = new byte[number];
for(int i = 0;i < number;i++)
{
    pattern[i] = z; // pattern you want to switch this channel to
}
device.SetChannels(pattern);</pre>
```

Get all channel switches at once:

```
byte[] pattern = device.GetChannels();
```

Set one channel switch:

```
ushort index = 10;
byte pattern = 1;
device.SetChannel(index, pattern)
```

Get one channel switch:

```
ushort index = 10;
byte pattern = device.GetChannel(index);
```

3 Function Classes 4

3 Function Classes

- Mcs.Usb.CCMOSMea_FunctionNet
- Mcs.Usb.CDacCalibrationFunctionNet
- Mcs.Usb.CDigOutStimulatorFunctionNet
- Mcs.Usb.CIntanMea_FunctionNet
- Mcs.Usb.CInterfaceboardFunctionNet
- Mcs.Usb.CMcsBus_MotorControlNet
- Mcs.Usb.CMcsBus_VoltageModeNet
- Mcs.Usb.CMcsBus_AxisParametersNet
- Mcs.Usb.CMcsBus SensorNet
- Mcs.Usb.CMcsBus_TempSensorNet
- Mcs.Usb.CMcsBus_ExtensionNet
- Mcs.Usb.CMcsBus_FYIExtensionNet
- Mcs.Usb.CMcsUsbDeviceStatePushFunctionNet
- Mcs.Usb.CMEA2100x256FunctionNet
- Mcs.Usb.CMeaAudioFunctionNet
- · Mcs.Usb.CMeaDigitalDataFunctionNet
- Mcs.Usb.CMeaFeedbackFunctionNet
- Mcs.Usb.CMeFunctionNet
- · Mcs.Usb.CMultiwellCallbackFunctionNet
- Mcs.Usb.CMultiwellOptoStimFunctionNet
- Mcs.Usb.CPPCFunctionNet
- Mcs.Usb.CPPS_FunctionNet
- Mcs.Usb.CPPS_FunctionNet
- Mcs.Usb.CPulseGeneratorFunctionNet
- Mcs.Usb.CRFFunctionNet
- Mcs.Usb.CRobo_FYITemp_FunctionNet
- Mcs.Usb.CRobo_FYIProgram_FunctionNet
- Mcs.Usb.CRobo_FYITemp_FunctionNet
- Mcs.Usb.CRobo_FYIProgram_FunctionNet
- Mcs.Usb.CSCUFunctionNet
- · Mcs.Usb.CStimulusFunctionNet
- Mcs.Usb.CTEERFunctionNet
- Mcs.Usb.CW2100_FunctionNet
- Mcs.Usb.CW2100_StimulatorFunctionNet

- Mcs.Usb.CWarnerUssingFunctionNet
- Mcs.Usb.CWarnerValveControllerDeviceTesterFunctionNet
- · Mcs.Usb.CWClassicFunctionNet
- · Mcs.Usb.CWirelessBaseFunctionNet

4 Data ACQuisition (DACQ) Devices

There are different device types of (MEA) data acquisition (DACQ) devices. All of them are supported by this class.

This library does **not** support the writing of the MCD (MC_Rack), MSRD (Multi Channel Experimenter) or HDF5 file format!

The class Mcs.Usb.CMeaDeviceNet is the base class for DACQ devices.

The base class Mcs.Usb.CMeaDeviceNet constructs actually the underlying classes for USB-MEA devices (Mcs.Usb.CMeaUSBDeviceNet).

```
CMeaDeviceNet device = new CMeaDeviceNet(McsBusTypeEnumNet.MCS_USB_BUS, OnChannelData, OnError);
```

For connecting to a DACQ device see Connecting to an MCS device.

Get the number of available analog hardware channels and set the number of channels to the maximum.

```
int hwchannels;
device.HWInfo().GetNumberOfHWADCChannels(out hwchannels);
device.SetNumberOfChannels(hwchannels);
int samplingrate = 1000;
device.SetSamplerate(samplingrate, 1, 0);
device.EnableDigitalIn(true, 0);
```

Get the layout to know how the data look like that you receive

```
int ana, digi, che, tim, block;
device.GetChannelLayout(out ana, out digi, out che, out tim, out block);
```

For the Mcs.Usb.OnChannelData callback function you have to provide a definition of the channels you want to receive.

```
bool[] selChannels = new bool[block];
for (int i = 0; i < block; i++)
{
    selChannels[i] = true; // With true channel i is selected
    // selChannels[i] = false; // With false the channel i is deselected
}
channelblocksize = samplingrate / 10;
// queue size and threshold should be selected carefully
device.SetSelectedChannels(selChannels, 10 * channelblocksize, channelblocksize);</pre>
```

The Mcs.Usb.OnChannelData callback function gets a callback for each channelblock that is defined. In this example a callback for each channel.

```
void OnChannelData(CMcsUsbDacqNet d, int cbHandle, int numSamples)
{
   int size_ret;
   ushort[] channeldata = device.ChannelBlock_ReadFramesUI16(CbHandle, numSamples, out size_ret);
}
void OnError(String msg, int info)
{
   MessageBox.Show("Mea Device Error: " + msg);
}
```

see MEA_Recording in the Examples directory.

5 The MCS Robo Device 6

5 The MCS Robo Device

5.1 Introduction

Up to now two MCS devices exist that base on the Robo platform.

- The MCS Roboinject device is controlled by the Mcs.Usb.CRobolnjectDeviceNet class.
- The MCS Roboocyte2 device is controlled by the Mcs.Usb.CRoboocyte2DeviceNet class.

Both classes are derived from Mcs.Usb.CRoboDeviceNet

6 STG200x & STG400x STimulus Generator

6.1 Introduction

The STG200x & STG400x Series Stimulus Generators have two distinct modes of operation, the Download mode and the Streaming mode.

6.2 Download mode

The Download mode is the "classic" mode of operation, as used by the MC Stimulus software. In this mode, one or multiple waveforms are defined in PC memory and downloaded to the STG. The waveforms are stored in STG device onboard memory and can be sent to the analog and sync outputs once or multiple times. The STG can operate independently from the PC (without computer connection) after the download. Output is triggered either by the front panel start/stop button, the digital trigger inputs or under software control.

In the Download mode, there are up to eight independent triggers available (depending on the device). The user can assign each of the analog outputs and sync (digital) outputs to any of the triggers.

The analog output waveform is stored sample by sample in the STG memory. To reduce memory usage, this data can be compressed: whenever a given output value is to be held for more than one sample period, it has only to be given once. The user can define the number of sample periods for that a pattern should remain active. Compression is done for each channel independently of the others, thus the algorithm to compress the data is very easy to implement.

A new feature of the Download mode is the segmentation of the STG memory. The onboard memory can be devided into up to 100 segments. Each segment can hold its own waveform pattern. Under software control, the user can switch between the defined segments within milliseconds. Another option is to use the four trigger inputs to select between four predefined segments. This option is accessible from the MC_Stimulus Software as the "Multi-File mode", and can start each of up to four defined waveforms within microseconds. This feature allows a predefinied flexible response (feedback) to recorded data.

Mcs.Usb.CStg200xDownloadNet is the class for using the STG in download mode.

6.2 Download mode 7

6.2.1 Memory Layout and Trigger Setup

The class to be used for the Download mode is Mcs.Usb.CStg200xDownloadNet, which is derived from Mcs.Usb.CStg200xBasicNet. You can add a poll handler delegate (Mcs.Usb.OnStg200xPollStatus) to the constructor Mcs.Usb.CStg200xDownloadNet.

For connecting to an STG see Connecting to an MCS device.

To use the Download mode, the memory layout of the STG200x can be set up, if the default is not sufficient. The total amount of memory available in the STG is obtained by the Mcs.Usb.CStg200xDownloadNet.GetTotalMemory call. With Mcs.Usb.CStg200xDownloadNet.SendSegmentDefine the segment sizes are assigned.

```
uint32_t memory = device.GetTotalMemory(); // obtain total memory available
uint[] segmentmemory = new uint[2]; // each segments has half of total memory
segmentmemory[0] = memory / 2;
segmentmemory[1] = memory / 2;
device.SendSegmentDefine(segmentmemory);// setup the STG
```

Next, for each segment, one has to assign the amount of memory to be used for each channel and sync output. This is done by Mcs.Usb.CStg200xDownloadBasicNet.SetCapacity. Its arguments contain a list of memory sizes, with one entry per channel and one entry per sync output. Again, the total memory assigned to the channels and sync outputs must not exceed the memory assigned to the segment.

Before the STG can start, the trigger has to be configured. This is done by the Mcs.Usb.CStg200xDownloadNet.SetupTrigger call. Its arguments are a list of channelmaps, syncoutmaps and repeats, one for each of the four available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and syncout 1 to trigger 1 and channel 3 to trigger 2 use:

```
uint32_t TriggerInputs = device.GetNumberOfTriggerInputs();
uint[] channelmap = new uint[TriggerInputs];
uint[] syncoutmap = new uint[TriggerInputs];
uint[] repeat = new uint[TriggerInputs];
for (int i = 0; i < TriggerInputs; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    repeat[i] = 0;
}
// Trigger 0
channelmap[0] = 1; // Channel 1
syncoutmap[0] = 1; // Syncout 1
repeat[0] = 0; // forever
// Trigger 1
channelmap[1] = 4; // Channel 3
device.SetupTrigger(channelmap, syncoutmap, repeat);</pre>
```

For the STG400x series you have to set the output mode of the channels. Mcs.Usb.CStg200xDownloadNet.SetVoltageMode interprets the values as voltages. Mcs.Usb.CStg200xDownloadNet.SetCurrentMode as currents.

```
// Only meaningfull for STG400x
device.SetVoltageMode();
```

For each segment, data can be sent to each of the defined channels and sync outputs using the Mcs.Usb.CStg200xDownloadNet.SendChannelData and Mcs.Usb.CStg200xDownloadNet.SendSyncData calls. channeldata and syncdata are a list of analog and digital samples as a list of two byte values (unsigned short). Multiple calls to Mcs.Usb.CStg200xDownloadNet.SendSyncData to the same channel append data to that channel.

6.3 Streaming mode 8

If the Multi-File mode of the STG is enabled using the Mcs.Usb.CStg200xDownloadNet.EnableMultiFileMode call, the four trigger inputs are used to switch between four segments. A hardware trigger signal (TTL) on trigger input 1 selects the first segment and starts all pulses in this segment. Thus with the Multi-File mode, one can predefine four stimulus patterns and switch between them without a connection to the PC.

The STG200x series has an analog resolution of 13 bits, thus the analog data contains the information in bits 0 to 12 of each sample. Bits 13 to 15 have to be 0.

```
int DACResolution = device.GetDACResolution();
// Data for Channel 0
   device.ClearChannelData(0);
   double factor = 0.1;
   const int 1 = 1000;
       ushort[] pData = new ushort[1];
       Uint64_t[] tData = new Uint64_t[1];
       for (int i = 0; i < 1; i++)
           // calculate Sin-Wave
           double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
               Math.Sin(2.0 * (double)i * Math.PI / (double)1);
           // calculate sign
           pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
               (int)Math.Pow(2, DACResolution - 1));
           tData[i] = (Uint64_t)20; // duration in \mus
       device.SendChannelData(0, pData, tData);
// Data for Channel 3
   device.ClearChannelData(2);
   double factor = 0.1;
   const int 1 = 700;
   // without compression
   ushort[] pData = new ushort[1];
   Uint64_t[] tData = new Uint64_t[1];
for (int i = 0; i < 1; i++)</pre>
       // calculate Sin-Wave
       double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
           Math.Sin(2.0 * (double)i * Math.PI / (double)1);
       // calculate sign
       device.SendChannelData(2, pData, tData);
// Data for Sync 0
   device.ClearSyncData(0);
   ushort[] pData = new ushort[1000];
   Uint64_t[] tData = new Uint64_t[1000];
   for (int i = 0; i < 1000; i++)
       pData[i] = (ushort)(i&1);
       tData[i] = 20;
   device.SendSyncData(0, pData, tData);
}
```

Start the trigger by pushing the front button or by software

```
// Start Trigger 1 and 2 device.SendStart(1 + 2); // Trigger 1 und 2 \,
```

see the StgDownloadExampleNet in the example directory.

6.3 Streaming mode

The other mode of operation is the Streaming mode. Here the analog output is sent to the STG device in "real time". The PC has to be connected to the STG all the time. The data that is sent to the analog output is downloaded from the PC to the STG on the fly.

The Streaming mode is useful for applications where flexible feedback is needed as well for applications where very long waveforms which are not repeated (such as white noise) are used.

6.3 Streaming mode 9

The Streaming mode works by use of two ring buffers which hold data. One is in PC memory and managed by the DLL, and one is in on-board STG memory. Data is transferred from PC memory to the STG via the USB bus in time slices of one millisecond.

The user can define both the size of the ring buffer in DLL memory and in the STG memory. Once the Streaming mode is started, the STG request data from the PC. The data rate from PC to STG is variable and controlled by the STG. The STG request data from the PC at a rate to keep its internal ringbuffer at about half full.

It is the responsibility of the user to keep the ring buffer in the memory of the PC filled, so the DLL can supply sufficient data to the STG. To do so, the Windows DLL allows to define a "callback" function which is called whenever new data is needed, or more precise, as soon as the ring buffer in the memory of the PC falls below the user defined threshold.

Small buffers have the advantage of a low latency between data generation in the callback funtion and its output as a analog signal from the STG. However for low latency to work, the user-written callback function has to be fast and to produce a steady flow of data.

In the Streaming mode, all triggers are available as well. Each of the eight analog and sync outputs can be assigned to one of the triggers.

The output rate is user defined with a maximum of 50 kHz

Mcs.Usb.CStg200xStreamingNet is the class for using the STG in streaming mode.

6.3.1 Memory Layout and Trigger Setup

With the constructor for Mcs.Usb.CStg200xStreamingNet.CStg200xStreamingNet, the name of the callback function for the data handler is provided. The data handler function is called automatically, whenever the STG needs new data. This data is first written to a ring buffer in the memory of the PC. The size for this ring buffer is defined as first argument in the constructor. The user provided delegate gets the trigger number which needs new data as argument

```
CStg200xStreamingNet device = new CStg200xStreamingNet(10000, dataHandler, errorHandler);
```

The callback funtion, which is defined in the constructor, is called whenever the STG needs new data for a trigger, or more precise, whenever the ring buffer in PC memory falls below the defined threshold.

The user can query the amount of space available for queuing by use of the Mcs.Usb.CStg200xStreamingNet.← GetDataQueueSpace call. Its return value is the number of samples that can be send to the STG.

User code is required to fill an array analog and sync out data, sample by sample for up to the maximum number of samples as obtained by Mcs.Usb.CStg200xStreamingNet.GetDataQueueSpace or Mcs.Usb.CStg200xStreaming← Net.GetSyncoutQueueSpace.

The values for the analog outputs are 16 bits signed integers. The lower bits are trunctated according to the resolution of the STG. This behaviour is different to the behaviour in download mode.

Note: Compression as described in the download mode can NOT be used for the streaming mode.

The new data is sent to the STG by using the Mcs.Usb.CStg200xStreamingNet.EngueueData call.

6.3 Streaming mode 10

```
for (int i = 0; i < 1000; i++)
                      // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
                      double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) *
    Math.Sin(2.0 * (double)i * Math.PI / (double)1000);
                      data[i] = (short)sin;
                 uint32_t enqueued = device.EnqueueData(channel, data);
         {// Handle Channel 3
             uint32_t channel = 2;
             for (; ; )
                 uint32_t space = device.GetDataQueueSpace(channel);
                 if (space < 700)</pre>
                 short[] data = new short[700];
                 for (int i = 0; i < 700; i++)
                      // Calc Sin-Wave (16 bits) lower bits will be removed according resolution
                      double sin = factor * (Math.Pow(2, 16 - 1) - 1.0)
                          Math.Sin(2.0 * (double)i * Math.PI / (double)700);
                      data[i] = (short)sin;
                 uint32_t enqueued = device.EnqueueData(channel, data);
         {// Handle Syncout 1
             uint32_t channel = 0;
             for (; ; )
                 uint32_t space = device.GetSyncoutQueueSpace(channel);
                 if (space < 1000)
                 ushort[] data = new ushort[1000];
                 for (int i = 0; i < 1000; i++)
  data[i] = (ushort)(i & 1);</pre>
                 uint32_t enqueued = device.EnqueueSyncout(channel, data);
        }
    }
void errorHandler()
```

For connecting to an STG device see Connecting to an MCS device.

With enabling or disabling the continuous mode it can be selected how the STG handles an "out of data" situation.

When Mcs.Usb.CStg200xStreamingNet.EnableContinousMode is used, the STG does not stop when it runs out of data, but it keeps running and sends a zero voltage to its outputs.

When Mcs.Usb.CStg200xStreamingNet.DisableContinousMode is used, the STG stops when it runs out of data. It has to be retriggered to resume the output.

device.EnableContinousMode();

Mcs.Usb.CStg200xStreamingNet.SetOutputRate is used to set the sampling rate. device.SetOutputRate(50000);

To use the Streaming mode, the memory layout of the STG has to be set up. To total amount of memory available in the STG is obtained by the Mcs.Usb.CStg200xStreamingNet.GetTotalMemory call.

This memory can be assigned to four ring buffers (one per trigger) which buffer the data received from the PC via USB cable. This is done with the CStg200xStreaming::SetCapacity call. The total amount of memory must not exceed the total memory size as obtained by Mcs.Usb.CStg200xStreamingNet.GetTotalMemory.

This internal ring buffer is crucial for proper operation of the Streaming mode. The size of the ring buffer determines the latency of the Streaming mode. The firmware of the STG requests data from the PC in order to keep the ring buffer about half full. Thus the average latency is:

```
latency = (ringbuffersize in bytes/4) / output rate
```

7 Namespace Index 11

If the ring buffer size is too big, the latency of the STG might be too long. If the ring buffer size is too low, an overflow or underflow of data in the STG ringbuffer might occur, resulting in data jumps of the output signals or the "out of data" situation described erlier.

The following example divides the total memory equally amoung the four triggers:

Before the STG can start, the trigger has to be configured. This is done by the Mcs.Usb.CStg200xStreaming
Net.SetupTrigger call. Its arguments are a list of channelmaps, syncoutmaps, digoutmap, autostart and callback
_threshold, with one entry for each of the available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and 3 and syncout 1 to trigger 1 use:

```
uint32_t ntrigger = device.GetNumberOfTriggerInputs();  // obtain number of triggers in this STG
uint[] channelmap = new uint[ntrigger];
uint[] syncoutmap = new uint[ntrigger];
uint[] digoutmap = new uint[ntrigger];
uint[] autostart = new uint[ntrigger];
uint[] callback_threshold = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    autostart[i] = 0;
    callback_threshold[i] = 0;
}
channelmap[0] = 0x1 + 0x4; // Channel 1 und Channel 3 to Trigger 1
syncoutmap[0] = 0x1; // Syncout 1 to Trigger 1
autostart[0] = 1;
callback_threshold[0] = 50; // 50% of buffer size
device.SetupTrigger(channelmap, syncoutmap, digoutmap, autostart, callback_threshold);
device.StartLoop();
System.Threading.Thread.Sleep(1000); // Give StartLoop some time</pre>
```

Start Trigger by pushing the front button or by Software

device.SendStart(1);

see the StgStreamingExampleNet in the example directory.

7 Namespace Index

7.1 Namespace List

Here is a list of all namespaces with brief descriptions:

Mcs::Usb 22

8 Hierarchical Index

8.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

8.1 Class Hierarchy 12

CW2100_FunctionNet::AudioChannelsNet	29
BatteryState	29
CCreateFilterNet	48
BesselFilterHighPassNet	30
BesselFilterLowPassNet	31
ButterworthFilterHighPassNet	31
ButterworthFilterLowPassNet	32
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf DacqGroupChannelEnumTemplateNet} >$	57
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf DacqGroupChannelEnumNet} >$	57
CDeviceGroupChannelInfoNet	56
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf int} >$	57
CDeviceGroupChannelInfoGenericNet	56
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf SCUDacqGroupChannelEnumNet} >$	57
CDeviceGroupChannelInfoSCUNet	57
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf W2100DacqGroupChannelEnumNet} >$	57
CDeviceGroupChannelInfoW2100Net	58
CFilterCoefficientsNet	66
CFilterPropertyNet	72
CMcsUsbDacqNet::CHWInfo	112
CMcsUsbFunctionNet	230
$\label{lem:condition} \textbf{CDacqGroupChannelSelectionTemplateNet} < \textbf{DacqGroupChannelEnumNet}, \textbf{DacqGroup} \leftarrow \\ \textbf{ChannelEnum, CDeviceGroupChannelInfoNet} >$	53
CDacqGroupChannelSelectionNet	52
${\tt CDacqGroupChannelSelectionTemplateNet} {< } {\tt int, int, CDeviceGroupChannelInfoGenericNet} {< }$	53
CDacqGroupChannelGenericSelectionNet	52
$\textbf{CDacqGroupChannelSelectionTemplateNet} < \textbf{SCUDacqGroupChannelEnumNet}, \textbf{SCUDacq} \leftarrow \textbf{GroupChannelEnum}, \textbf{CDeviceGroupChannelInfoSCUNet} >$	53
CSCUDacqGroupChannelSelectionNet	420
$\label{lem:condition} \textbf{CDacqGroupChannelSelectionTemplateNet} < & \textbf{W2100DacqGroupChannelEnumNet}, & \textbf{W2100} \leftarrow \textbf{DacqGroupChannelEnum}, & \textbf{CDeviceGroupChannelInfoW2100Net} > \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ \\ $	53
CW2100DacqGroupChannelSelectionNet	539
CCMOSMea_FunctionNet	34
CDacCalibrationFunctionNet	50

CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, GroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >	Dacq <i>←</i> 53
CDigOutStimulatorFunctionNet	59
CFilterConfigurationNet	68
CFilterConfigurationRegisterNet	70
CGrapheneFunctionNet	102
CIntanMea_FunctionNet	114
CInterfaceboardFunctionNet	118
CInterfaceboard2FunctionNet	116
CMcsBus_AxisParametersNet	129
CMcsBus_ExtensionNet	130
CMcsBus_FYIExtensionNet	131
CMcsBus_MotorControlNet	133
CMcsBus_SensorNet	150
CMcsBus_TempSensorNet	161
CMcsBus_VoltageModeNet	163
CMcsBusNet	167
CMcsUsbDeviceStatePushFunctionNet	219
CMultiwellCallbackFunctionNet	304
CSCUFunctionNet	421
CMEA2100x256FunctionNet	261
CMeaAudioFunctionNet	263
CMeaDigitalDataFunctionNet	284
CMeaFeedbackFunctionNet	286
CMeFunctionNet	295
CMultiwellOptoStimFunctionNet	312
CPPCFunctionNet	345
CPPS_FunctionNet	355
CProgramPressureCurveNet	360
CPulseGeneratorFunctionNet	361
CRFFunctionNet	369
CRobo FYIProgram FunctionNet	374

CRobo_FYITemp_FunctionNet	375
CStimulusFunctionNet	484
CTEERFunctionNet	511
CW2100_StimulatorFunctionNet	531
CWarnerUssingFunctionNet	541
CWarnerValveControllerDeviceTesterFunctionNet	580
CWirelessBaseFunctionNet	588
CW2100_FunctionNet	524
CWClassicFunctionNet	583
CMcsUsbFunctionPointerContainer	231
CMcsUsbListEntryNet	231
CMcsUsbListNet	237
CMcsUsbNet	240
CExternDTesterDeviceNet	64
CFluidControlDeviceNet	73
CGenericDevelopDeviceNet	82
CGilsonDeviceNet	100
CMcsUsbDacqNet	171
CMeaDeviceNet	276
CMeaUSBDeviceNet	294
CCMOSMeaDeviceNet	45
CHLADacqNet	110
CLIH3DeviceNet	120
CMultiwellDeviceNet	306
CWarnerUssingDeviceNet	540
COctoPotDeviceNet	318
CRoboDacqNet	377
CMcsUsbDeviceStatePushNet	220
CWarnerValveControllerDeviceNet	556
CMcsUsbFactoryNet	221
CMeaCleanDeviceNet	266
CMeaCoatDeviceNet	270

CMealmpedanceDeviceNet	290
CMeaSwitchDeviceNet	292
CChannelTestDeviceNet	33
CMultiBatteryChargerDeviceNet	297
CNF_GenDeviceNet	317
COkuvisionStimulatorDeviceNet	322
CPathIdentDeviceNet	326
CPedoterDeviceNet	327
CPeristalticPumpDeviceNet	329
CPgaDeviceNet	330
CPositionIIDeviceNet	332
CPositionImpDeviceNet	341
CPPCDeviceNet	344
CPPS_DeviceNet	354
CRadioControledDevicesNet	364
CRetinaLedDeviceNet	367
CRoboDeviceNet	391
CEncapsulatorDeviceNet	63
CFYIDeviceNet	80
CHiClampDeviceNet	109
CHLADeviceNet	111
CMeasureTableDeviceNet	291
CPatchServerDeviceNet	325
CPPSDeviceNet	359
CRoboStatorDeviceNet	411
CRobolnjectDeviceNet	409
CRoboocyte2DeviceNet	409
CTEERMachineDeviceNet	521
CRoboFluidDeviceNet	405
CSafeISDeviceNet	417
CSerialPortNet	436
CStg200xBasicNet	438

9 Class Index 16

CStg200xDownloadBasicNet	468
CStg200xDownloadNet	477
CSw2to64DeviceNet	495
CTcxDeviceNet	498
CMcsUsbPointerContainer	261
CCMOSMeaDeviceNet::CRegionOfInterestRect	366
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	523
DeviceIdNet	588
DigitalSource< digitalsourceenum >	590
DriverVersionNet Exception	591
CUsbExceptionNet	522
FirmwareDestinationNames	596
HeadstageIDTypeObject	603
HeadStageIDTypeState IComparable	604
HeadStageIDType	600
mkfilterNet	605
CRoboDeviceNet::RoboMainLowLevelCommands	608
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	615
CFilterCoefficientsNet::s_FilterAttributesNet	615
CMeaAudioFunctionNet::s_setaudionet	616
CStimulusFunctionNet::SidebandData	617
StgStatusNet stgstreaming	618
CStg200xBasicNet	438
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	618
usbSetupPacket_t	620
W2100_StimulusParametersNet	620

9 Class Index

9.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

CW2100_FunctionNet::AudioChannelsNet	29
BatteryState	29
BesselFilterHighPassNet	30
BesselFilterLowPassNet	31
ButterworthFilterHighPassNet	31
ButterworthFilterLowPassNet	32
CChannelTestDeviceNet	33
CCMOSMea_FunctionNet	34
CCMOSMeaDeviceNet	45
CCreateFilterNet	48
CDacCalibrationFunctionNet 50	
CDacqGroupChannelGenericSelectionNet	52
CDacqGroupChannelSelectionNet	52
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupCha 53	nnelEnumTemplate
CDeviceGroupChannelInfoGenericNet	56
CDeviceGroupChannelInfoNet	56
CDeviceGroupChannelInfoSCUNet	57
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf DacqGroupChannelEnumTemplateNet} >$	57
CDeviceGroupChannelInfoW2100Net	58
CDigOutStimulatorFunctionNet CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class	59
CEncapsulatorDeviceNet CEncapsulatorDeviceNet is the to control the MCS HiClamp device	63
CExternDTesterDeviceNet CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)	64
CFilterCoefficientsNet	66
CFilterConfigurationNet	68
CFilterConfigurationRegisterNet	70
CFilterPropertyNet	72
CFluidControlDeviceNet CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device	73
CFYIDeviceNet CFYIDeviceNet is the class to control the MCS FYI device	80

CGenericDevelopDeviceNet CGenericDevelopDeviceNet is the class to use during development of a new device	82
CGilsonDeviceNet	
CGilsonDeviceNet is the class to control a Gilson device	100
CGrapheneFunctionNet CGrapheneFunctionNet is the class to control the TEER device	102
CHiClampDeviceNet CHiClampDeviceNet is the to control the MCS HiClamp device	109
CHLADacqNet	110
CHLADeviceNet CHLADeviceNet is the to control the MCS HLA device	111
CMcsUsbDacqNet::CHWInfo Class to provide hardware information about the device	112
CIntanMea_FunctionNet	114
CInterfaceboard2FunctionNet CInterfaceboard2FunctionNet is the class to control the Interfaceboard	116
CInterfaceboardFunctionNet CInterfaceboardFunctionNet is the class to control the Interfaceboard	118
CLIH3DeviceNet CLIH3DeviceNet is the class to access the HEKA LIH3 device	120
CMcsBus_AxisParametersNet	129
CMcsBus_ExtensionNet	130
CMcsBus_FYIExtensionNet	131
CMcsBus_MotorControlNet	133
CMcsBus_SensorNet	150
CMcsBus_TempSensorNet	161
CMcsBus_VoltageModeNet	163
CMcsBusNet	167
CMcsUsbDacqNet Base class for data acquisition devices	171
CMcsUsbDeviceStatePushFunctionNet	219
CMcsUsbDeviceStatePushNet	220
CMcsUsbFactoryNet	221
CMcsUsbFunctionNet	230
CMcsUsbFunctionPointerContainer	231
CMcsUsbListEntryNet McsUsbListEntryNet identifies a connected device	231

CMcsUsbListNet Class to handle a list of connected MCS USB devices	237
CMcsUsbNet Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class	240
CMcsUsbPointerContainer	26 1
CMEA2100x256FunctionNet CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include " Stg200xNet.h" to resolve documentation reference	26 1
CMeaAudioFunctionNet	263
CMeaCleanDeviceNet CMeaCleanDeviceNet is the class to access the MEA Clean device	266
CMeaCoatDeviceNet CMeaCoatDeviceNet is the class to access the MEA Coat device	270
CMeaDeviceNet Base class for MEA data acquisition devices	276
CMeaDigitalDataFunctionNet	284
CMeaFeedbackFunctionNet	286
CMealmpedanceDeviceNet	290
CMeasureTableDeviceNet CMeasureTableDeviceNet is the to control the MCS HLA device	29 1
CMeaSwitchDeviceNet The class to control the USB-MEA-Switch	292
CMeaUSBDeviceNet Class for data acquisition via ME and MEA USB amplifiers	294
CMeFunctionNet 295	
CMultiBatteryChargerDeviceNet CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device	297
CMultiwellCallbackFunctionNet CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator	304
CMultiwellDeviceNet CMultiwellDeviceNet is the class to access the Multiwell device	306
CMultiwellOptoStimFunctionNet CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device	312
CNF_GenDeviceNet	317
COctoPotDeviceNet	318
COkuvisionStimulatorDeviceNet	322

CPatchServerDeviceNet CPatchServerDeviceNet is the class to control the MCS PatchServer device	325
CPathIdentDeviceNet CPathIdentDeviceNet	326
CPedoterDeviceNet 327	
CPeristalticPumpDeviceNet CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump	329
CPgaDeviceNet	330
CPositionIIDeviceNet CPositionIIDeviceNet is the class to control PositionII devices	332
CPositionImpDeviceNet CPositionImpDeviceNet is the class to access the Position/Imp devices	341
CPPCDeviceNet	344
CPPCFunctionNet CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump	345
CPPS_DeviceNet	354
CPPS_FunctionNet	355
CPPSDeviceNet CPPS4plus1DeviceNet is the to control the MCS HLA device	359
CProgramPressureCurveNet is the class to program pressure curves	360
CPulseGeneratorFunctionNet CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking	361
CRadioControledDevicesNet	364
CCMOSMeaDeviceNet::CRegionOfInterestRect	366
CRetinaLedDeviceNet	367
CRFFunctionNet CRFFunctionNet is the class to control RF devices	369
CRobo_FYIProgram_FunctionNet	374
CRobo_FYITemp_FunctionNet	375
CRoboDacqNet	377
CRoboDeviceNet CRoboDeviceNet is the base class for all Robo platform based devices	391
CRoboFluidDeviceNet	405
CRobolnjectDeviceNet CRobolnjectDeviceNet is the to control the MCS Robolnject device	409
CRoboocyte2DeviceNet CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device	409

CRoboStatorDeviceNet	411
CSafelSDeviceNet 417	
CSCUDacqGroupChannelSelectionNet	420
CSCUFunctionNet CSCUFunctionNet is the class to control the SCU device	421
CSerialPortNet	436
CStg200xBasicNet Base class for the Stg200x	438
CStg200xDownloadBasicNet CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device	468
CStg200xDownloadNet Main class for the STG download mode This class implements the STG download mode interface.	477
CStimulusFunctionNet	484
CSw2to64DeviceNet The class to control the MCS-USB-Sw2to64 device	495
CTcxDeviceNet Class to control a Temperature Controller (TCX)	498
CTEERFunctionNet CTEERFunctionNet is the class to control the TEER device	511
CTEERMachineDeviceNet	52 1
CUsbExceptionNet Exception class that is thrown in case of an USB error	522
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	523
CW2100_FunctionNet	52 4
CW2100_StimulatorFunctionNet	53 1
CW2100DacqGroupChannelSelectionNet	539
CWarnerUssingDeviceNet CWarnerUssingDeviceNet is the class to control the Ussing device	540
CWarnerUssingFunctionNet CWarnerUssingFunctionNet is the class to control the Ussing device	541
CWarnerValveControllerDeviceNet CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller	556
CWarnerValveControllerDeviceTesterFunctionNet CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester	580
CWClassicFunctionNet	583

CWirelessBaseFunctionNet	588
DeviceIdNet Device Id	588
DigitalSource< digitalsourceenum >	590
DriverVersionNet Class gives firmware versions of the device's firmware destinations	591
FirmwareDestinationNames	596
HeadStageIDType	600
HeadstageIDTypeObject	603
HeadStageIDTypeState	604
mkfilterNet	605
CRoboDeviceNet::RoboMainLowLevelCommands	608
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	615
CFilterCoefficientsNet::s_FilterAttributesNet	615
CMeaAudioFunctionNet::s_setaudionet	616
CStimulusFunctionNet::SidebandData	617
StgStatusNet	618
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	618
usbSetupPacket_t	620
W2100_StimulusParametersNet	620

10 Namespace Documentation

10.1 Mcs Namespace Reference

Namespaces

• Usb

10.2 Mcs::Usb Namespace Reference

Classes

- class BatteryState
- class BesselFilterHighPassNet
- class BesselFilterLowPassNet
- class ButterworthFilterHighPassNet
- class ButterworthFilterLowPassNet
- class CChannelTestDeviceNet

- class CCMOSMea FunctionNet
- class CCMOSMeaDeviceNet
- class CCreateFilterNet
- class CDacCalibrationFunctionNet
- class CDacqGroupChannelGenericSelectionNet
- class CDacqGroupChannelSelectionNet
- class CDacqGroupChannelSelectionTemplateNet
- class CDeviceGroupChannelInfoGenericNet
- class CDeviceGroupChannelInfoNet
- · class CDeviceGroupChannelInfoSCUNet
- class CDeviceGroupChannelInfoTemplateNet
- class CDeviceGroupChannelInfoW2100Net
- · class CDigOutStimulatorFunctionNet

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

· class CEncapsulatorDeviceNet

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

class CExternDTesterDeviceNet

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

- · class CFilterCoefficientsNet
- class CFilterConfigurationNet
- · class CFilterConfigurationRegisterNet
- class CFilterPropertyNet
- class CFluidControlDeviceNet

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

class CFYIDeviceNet

CFYIDeviceNet is the class to control the MCS FYI device

class CGenericDevelopDeviceNet

CGenericDevelopDeviceNet is the class to use during development of a new device.

• class CGilsonDeviceNet

CGilsonDeviceNet is the class to control a Gilson device.

class CGrapheneFunctionNet

CGrapheneFunctionNet is the class to control the TEER device

· class CHiClampDeviceNet

CHiClampDeviceNet is the to control the MCS HiClamp device

- class CHLADacqNet
- class CHLADeviceNet

CHLADeviceNet is the to control the MCS HLA device

- · class CIntanMea FunctionNet
- · class CInterfaceboard2FunctionNet

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

· class CInterfaceboardFunctionNet

CInterfaceboardFunctionNet is the class to control the Interfaceboard

· class CLIH3DeviceNet

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

- class CMcsBus_AxisParametersNet
- class CMcsBus_ExtensionNet
- class CMcsBus FYIExtensionNet
- · class CMcsBus MotorControlNet
- · class CMcsBus_SensorNet
- · class CMcsBus_TempSensorNet
- · class CMcsBus_VoltageModeNet
- class CMcsBusNet

· class CMcsUsbDacqNet

Base class for data acquisition devices.

- class CMcsUsbDeviceStatePushFunctionNet
- class CMcsUsbDeviceStatePushNet
- · class CMcsUsbFactoryNet
- class CMcsUsbFunctionNet
- · class CMcsUsbFunctionPointerContainer
- class CMcsUsbListEntryNet

McsUsbListEntryNet identifies a connected device.

class CMcsUsbListNet

Class to handle a list of connected MCS USB devices.

class CMcsUsbNet

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

- class CMcsUsbPointerContainer
- class CMEA2100x256FunctionNet

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

- class CMeaAudioFunctionNet
- class CMeaCleanDeviceNet

CMeaCleanDeviceNet is the class to access the MEA Clean device.

class CMeaCoatDeviceNet

CMeaCoatDeviceNet is the class to access the MEA Coat device.

class CMeaDeviceNet

Base class for MEA data acquisition devices.

- class CMeaDigitalDataFunctionNet
- · class CMeaFeedbackFunctionNet
- · class CMealmpedanceDeviceNet
- class CMeasureTableDeviceNet

CMeasure Table Device Net is the to control the MCS HLA device

• class CMeaSwitchDeviceNet

The class to control the USB-MEA-Switch.

· class CMeaUSBDeviceNet

Class for data acquisition via ME and MEA USB amplifiers

- class CMeFunctionNet
- · class CMultiBatteryChargerDeviceNet

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

· class CMultiwellCallbackFunctionNet

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

· class CMultiwellDeviceNet

CMultiwellDeviceNet is the class to access the Multiwell device.

class CMultiwellOptoStimFunctionNet

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

- · class CNF_GenDeviceNet
- · class COctoPotDeviceNet
- · class COkuvisionStimulatorDeviceNet
- · class CPatchServerDeviceNet

CPatchServerDeviceNet is the class to control the MCS PatchServer device

- class CPathIdentDeviceNet
- · class CPedoterDeviceNet
- class CPeristalticPumpDeviceNet

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

- · class CPgaDeviceNet
- · class CPositionIIDeviceNet

CPositionIIDeviceNet is the class to control PositionII devices

class CPositionImpDeviceNet

CPositionImpDeviceNet is the class to access the Position/Imp devices

- class CPPCDeviceNet
- · class CPPCFunctionNet

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

- class CPPS DeviceNet
- class CPPS FunctionNet
- · class CPPSDeviceNet

CPPS4plus1DeviceNet is the to control the MCS HLA device

· class CProgramPressureCurveNet

CProgramPressureCurveNet is the class to program pressure curves

class CPulseGeneratorFunctionNet

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

- · class CRadioControledDevicesNet
- · class CRetinaLedDeviceNet
- · class CRFFunctionNet

CRFFunctionNet is the class to control RF devices

- class CRobo_FYIProgram_FunctionNet
- class CRobo_FYITemp_FunctionNet
- class CRoboDacqNet
- · class CRoboDeviceNet

CRoboDeviceNet is the base class for all Robo platform based devices

- · class CRoboFluidDeviceNet
- class CRobolnjectDeviceNet

CRobolnjectDeviceNet is the to control the MCS Robolnject device

class CRoboocyte2DeviceNet

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

- · class CRoboStatorDeviceNet
- · class CSafeISDeviceNet
- class CSCUDacqGroupChannelSelectionNet
- · class CSCUFunctionNet

CSCUFunctionNet is the class to control the SCU device

- · class CSerialPortNet
- class CStg200xBasicNet

Base class for the Stg200x.

class CStg200xDownloadBasicNet

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

• class CStg200xDownloadNet

Main class for the STG download mode This class implements the STG download mode interface.

- class CStimulusFunctionNet
- · class CSw2to64DeviceNet

The class to control the MCS-USB-Sw2to64 device.

class CTcxDeviceNet

Class to control a Temperature Controller (TCX)

class CTEERFunctionNet

CTEERFunctionNet is the class to control the TEER device

- · class CTEERMachineDeviceNet
- class CUsbExceptionNet

Exception class that is thrown in case of an USB error.

- · class CW2100 FunctionNet
- class CW2100_StimulatorFunctionNet
- class CW2100DacqGroupChannelSelectionNet
- · class CWarnerUssingDeviceNet

CWarnerUssingDeviceNet is the class to control the Ussing device

class CWarnerUssingFunctionNet

CWarnerUssingFunctionNet is the class to control the Ussing device

class CWarnerValveControllerDeviceNet

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

class CWarnerValveControllerDeviceTesterFunctionNet

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

- class CWClassicFunctionNet
- class CWirelessBaseFunctionNet
- struct DeviceIdNet

Device Id.

- · class DigitalSource
- class DriverVersionNet

Class gives firmware versions of the device's firmware destinations.

- class FirmwareDestinationNames
- · class HeadStageIDType
- · class HeadstageIDTypeObject
- · class HeadStageIDTypeState
- class mkfilterNet
- class StgStatusNet
- class usbSetupPacket_t
- struct W2100_StimulusParametersNet

Enumerations

```
    enum enCMosMeaChipType {
        unknown = 0,
        nMos16LV = 1,
        nMos32LV = 3,
        nMos36LN = 6,
        nMos64LN = 7 }
    enum EnSTG200x_STATUS {
        OK,
        NOT_CONNECTED,
        DEVICE_NOT_FOUND }
```

Functions

- public delegate void OnMcsUsbDeviceState (usbSetupPacket t[^] request)
- private delegate void OnMcsUsbDeviceStateCallback (IntPtr pThis, uint32 t size, IntPtr buffer)
- public delegate void OnUpdateFirmwareStatusChange (String[^])
- public delegate void OnUpdateFirmwareProgress (int)
- public delegate void OnDeviceArrivalRemoval (CMcsUsbListEntryNet[^] entry)

Delegate to show a device arrival or removal.

- public delegate void OnStgPollStatus (unsigned int status, StgStatusNet[^] stgStatusNet, array< int >[^] index list)
- public delegate void OnMwPollStatus (unsigned int CurrentTemp, unsigned int PlateState, unsigned int SwitchState)

- public delegate void RoboStatusEventDelegate (array< unsigned char >^ buffer)
- public delegate void OnStg200xDataHandler (uint32_t trigger)
- public delegate void OnStg200xErrorHandler ()
- public delegate void OnChannelData (CMcsUsbDacqNet^ dacq, int CbHandle, int numFrames)
- public delegate void OnError (String[^] msg, int action)

10.2.1 Enumeration Type Documentation

10.2.1.1 enCMosMeaChipType enum enCMosMeaChipType [strong]

Enumerator

unknown	
nMos16LV	
nMos32LV	
nMos36LN	
nMos64LN	

10.2.1.2 EnSTG200x_STATUS enum EnSTG200x_STATUS [strong]

Enumerator

OK	
NOT_CONNECTED	
DEVICE_NOT_FOUND	

10.2.2 Function Documentation

```
10.2.2.2 OnDeviceArrivalRemoval() public delegate void Mcs::Usb::OnDeviceArrivalRemoval (
CMcsUsbListEntryNet^ entry)
```

Delegate to show a device arrival or removal.

```
10.2.2.3 OnError() public delegate void Mcs::Usb::OnError (
             String^ msg,
             int action )
10.2.2.4 OnMcsUsbDeviceState() public delegate void OnMcsUsbDeviceState (
             usbSetupPacket_t^ request )
10.2.2.5 OnMcsUsbDeviceStateCallback() private delegate void OnMcsUsbDeviceStateCallback (
             IntPtr pThis,
             uint32_t size,
             IntPtr buffer )
10.2.2.6 OnMwPollStatus() public delegate void Mcs::Usb::OnMwPollStatus (
             unsigned int CurrentTemp,
             unsigned int PlateState,
             unsigned int SwitchState )
10.2.2.7 OnStg200xDataHandler() public delegate void Mcs::Usb::OnStg200xDataHandler (
             uint32_t trigger )
10.2.2.8 OnStg200xErrorHandler() public delegate void Mcs::Usb::OnStg200xErrorHandler ( )
10.2.2.9 OnStgPollStatus() public delegate void Mcs::Usb::OnStgPollStatus (
             unsigned int status,
             StgStatusNet^ stgStatusNet,
             array < int >^{\wedge} index_list)
10.2.2.10 OnUpdateFirmwareProgress() public delegate void Mcs::Usb::OnUpdateFirmwareProgress
             int )
```

11 Class Documentation 29

10.2.2.11 OnUpdateFirmwareStatusChange() public delegate void Mcs::Usb::OnUpdateFirmware←

11 Class Documentation

11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference

Public Attributes

- W2100DacqGroupChannelEnumNet dacqgroup
- int channel
- · int amplification

11.1.1 Member Data Documentation

```
11.1.1.1 amplification int amplification
```

11.1.1.2 channel int channel

11.1.1.3 dacqgroup W2100DacqGroupChannelEnumNet dacqgroup

11.2 BatteryState Class Reference

Properties

```
unsigned int Charge [get]
unsigned int Voltage [get]
System::String^ ChargeString [get]
System::String^ ChargeRegionString [get]
```

• System::String^ VoltageString [get]

Generated by Doxygen

11.2.1 Property Documentation

11.2.1.1 Charge unsigned int Charge [get]

11.2.1.2 ChargeRegionString System:: String^ ChargeRegionString [get]


11.2.1.3 ChargeString System:: String ChargeString [get]

11.2.1.4 Voltage unsigned int Voltage [get]

11.2.1.5 VoltageString System:: String^ VoltageString [get]

11.3 BesselFilterHighPassNet Class Reference

Inheritance diagram for BesselFilterHighPassNet:

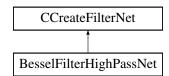

Public Member Functions

• BesselFilterHighPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.3.1 Constructor & Destructor Documentation

11.3.1.1 BesselFilterHighPassNet() BesselFilterHighPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.4 BesselFilterLowPassNet Class Reference

Inheritance diagram for BesselFilterLowPassNet:

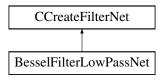

Public Member Functions

BesselFilterLowPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.4.1 Constructor & Destructor Documentation

11.4.1.1 BesselFilterLowPassNet() BesselFilterLowPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.5 ButterworthFilterHighPassNet Class Reference

Inheritance diagram for ButterworthFilterHighPassNet:

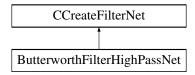

Public Member Functions

ButterworthFilterHighPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)


11.5.1 Constructor & Destructor Documentation

11.5.1.1 ButterworthFilterHighPassNet() ButterworthFilterHighPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.6 ButterworthFilterLowPassNet Class Reference

Inheritance diagram for ButterworthFilterLowPassNet:

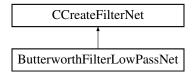

Public Member Functions

ButterworthFilterLowPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.6.1 Constructor & Destructor Documentation

11.6.1.1 ButterworthFilterLowPassNet() ButterworthFilterLowPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.7 CChannelTestDeviceNet Class Reference

Inheritance diagram for CChannelTestDeviceNet:

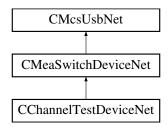

Public Member Functions

- CChannelTestDeviceNet ()
- ∼CChannelTestDeviceNet ()
- void SetWaveform (unsigned int Waveform)
- void SetAmplitude (unsigned int Amplitude)
- void SetFrequency (unsigned int Frequency)
- · void SetAttenuation (unsigned int Attenuation)

Additional Inherited Members

11.7.1 Constructor & Destructor Documentation


```
11.7.1.1 CChannelTestDeviceNet() CChannelTestDeviceNet ( )
```

```
11.7.1.2 \sim CChannelTestDeviceNet() \sim CChannelTestDeviceNet ( )
```

11.7.2 Member Function Documentation

11.8 CCMOSMea FunctionNet Class Reference

Inheritance diagram for CCMOSMea_FunctionNet:

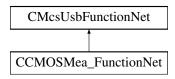

Public Member Functions

- CCMOSMea_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] cMOSMea_←
 FunctionPointerContainer)
- CCMOSMea_FunctionNet (CMcsUsbNet^ mcsusb)
- void SetADCInputOffset (int32_t offset)
- int32_t GetADCInputOffset ()
- void SetSourceDrain (int32_t voltage)
- int32 t GetSourceDrain ()
- void SetSourceGate (int32_t voltage)
- int32_t GetSourceGate ()
- void SetSourceBulk (int32_t voltage)
- int32 t GetSourceBulk ()
- void SetGate (int32_t voltage)
- int32 t GetGate ()
- void SetBath (int32_t voltage)
- int32_t GetBath ()
- int32 t GetGNDI ()
- int32_t GetVDDI ()
- int32_t GetVDD3I ()
- void UpdateTransistorVoltages ()
- bool AreTransistorVoltagesSet ()
- void PowerChip (bool on)
- bool IsChipPowered ()
- enCMosMeaChipType DetectChipType ()
- void SetGateToVOP ()
- void SetGateFloating ()
- bool IsGateFloating ()
- void VOPSTimerSetResetTimes (uint32_t ResetTime, uint32_t IntervalTime)
- void VOPSTimerSetResetTimes (uint32_t ResetTime, uint32_t IntervalTime, uint32_t HPFilterResetTime)
- void SetBathMode (CMOSMeaBathModeEnumNet Mode)
- CMOSMeaBathModeEnumNet GetBathMode ()
- void SetNeurochipMemoryData (uint16_t MemAddress, uint32_t MemData)

- void SetNeurochipMemoryData (uint16_t MemAddress, array< uint32_t >^ MemData)
- uint32_t GetNeurochipMemoryData (uint16_t MemAddress)
- array< uint32 t > ^ GetNeurochipMemoryData (uint16 t MemAddress, uint32 t RegestLength)
- uint32 t GetNeurochipMemorySize ()
- uint32_t GetMaxNumOfColumns (uint32_t Samplerate)
- void SetStimulusSites (List< int16_t >^ SwitchPosition)
- List< int16 t > ^ GetStimulusSites ()
- void ClearSTGOutput (uint32 t Channel)
- uint32 t GetNumberOfSupportedGroups ()
- uint32 t GetNumberOfSupportedGroups (uint32 t virtualDevice)
- DacqGroupChannelEnumNet GetGroupID (uint32 t Index)
- DacqGroupChannelEnumNet GetGroupID (uint32_t Index, uint32_t virtualDevice)
- uint32_t GetGroupNumberOfChannels (DacqGroupChannelEnumNet GroupID)
- uint32_t GetGroupNumberOfChannels (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- DacqMeaGroupTypeEnumNet GetGroupType (DacqGroupChannelEnumNet GroupID)
- void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBit← Map)
- void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >[^] EnabledChannelsBit←
 Map, uint32 t virtualDevice)
- List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumNet GroupID)
- List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- uint32 t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID)
- uint32 t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID)
- CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- int32_t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID)
- int32_t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- int32_t GetGroupADCBits (DacqGroupChannelEnumNet GroupID)
- int32 t GetGroupADCBits (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel← Number)
- uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel
 — Number, uint32_t virtualDevice)
- CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber)
- CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMealFDigChannelEnumNet GetGroupChannelBitmaskIFDigChannels (uint32 t ChannelNumber)
- CMOSMealFDigChannelEnumNet GetGroupChannelBitmasklFDigChannels (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel
 — Number)
- CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel
 — Number, uint32_t virtualDevice)

- CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t Channel↔ Number)
- CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t Channel
 — Number, uint32_t virtualDevice)
- CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32_t ChannelNumber)
- CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel
 — Number)
- CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel Number, uint32_t virtualDevice)
- CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_← t ChannelNumber)
- CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_← t ChannelNumber, uint32_t virtualDevice)

11.8.1 Constructor & Destructor Documentation

```
11.8.2.2 ClearSTGOutput() void ClearSTGOutput ( uint32_t Channel)
```

11.8.2.3 DetectChipType() enCMosMeaChipType DetectChipType ()

```
11.8.2.4 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (
              DacqGroupChannelEnumNet GroupID,
              List<br/>< bool >^ EnabledChannelsBitMap )
\textbf{11.8.2.5} \quad \textbf{EnableChannelsInGroup() [2/2]} \quad \texttt{void EnableChannelsInGroup ()}
              DacqGroupChannelEnumNet GroupID,
              List < bool >^{\wedge} EnabledChannelsBitMap,
              uint32_t virtualDevice )
11.8.2.6 GetADCInputOffset() int32_t GetADCInputOffset ( )
11.8.2.7 GetBath() int32_t GetBath ()
11.8.2.8 GetBathMode() CMOSMeaBathModeEnumNet GetBathMode ( )
11.8.2.9 GetEnabledChannelsInGroup() [1/2] List<br/>bool> ^ GetEnabledChannelsInGroup (
              DacqGroupChannelEnumNet GroupID )
11.8.2.10 GetEnabledChannelsInGroup() [2/2] List<bool> ^ GetEnabledChannelsInGroup (
              DacqGroupChannelEnumNet GroupID,
              uint32_t virtualDevice )
11.8.2.11 GetGate() int32_t GetGate ()
11.8.2.12 GetGNDI() int32_t GetGNDI ()
11.8.2.13 GetGroupADCBits() [1/2] int32_t GetGroupADCBits (
              {\tt DacqGroupChannelEnumNet} \  \, \textit{GroupID} \  \, )
```

```
11.8.2.14 GetGroupADCBits() [2/2] int32_t GetGroupADCBits (
            DacqGroupChannelEnumNet GroupID,
            uint32_t virtualDevice )
11.8.2.15 GetGroupChannelBitmaskBySelect() [1/2] uint32_t GetGroupChannelBitmaskBySelect (
            DacqGroupChannelEnumNet GroupID,
            uint32_t ChannelNumber )
11.8.2.16 GetGroupChannelBitmaskBySelect() [2/2] uint32_t GetGroupChannelBitmaskBySelect (
            DacqGroupChannelEnumNet GroupID,
            uint32_t ChannelNumber,
            uint32_t virtualDevice )
11.8.2.17 GetGroupChannelBitmaskHS1NCBathCurrent() [1/2] CMOSMeaHeadstage1NCBathCurrentEnum
Net GetGroupChannelBitmaskHS1NCBathCurrent (
            uint32_t ChannelNumber )
11.8.2.18 GetGroupChannelBitmaskHS1NCBathCurrent() [2/2] CMOSMeaHeadstagelNCBathCurrentEnum↔
Net GetGroupChannelBitmaskHS1NCBathCurrent (
            uint32_t ChannelNumber,
            uint32_t virtualDevice )
11.8.2.19 GetGroupChannelBitmaskHS1NCCol2Current() [1/2] CMOSMeaHeadstagelNCCol2CurrentEnum←
Net GetGroupChannelBitmaskHS1NCCol2Current (
            uint32_t ChannelNumber )
11.8.2.20 GetGroupChannelBitmaskHS1NCCol2Current() [2/2] CMOSMeaHeadstage1NCCol2CurrentEnum←
Net GetGroupChannelBitmaskHS1NCCol2Current (
            uint32_t ChannelNumber,
            uint32_t virtualDevice )
11.8.2.21 GetGroupChannelBitmaskHS1NChipTemp() [1/2] CMOSMeaHeadstage1NChipTempEnumNet Get↔
GroupChannelBitmaskHS1NChipTemp (
            uint32_t ChannelNumber )
```

```
11.8.2.22 GetGroupChannelBitmaskHS1NChipTemp() [2/2] CMOSMeaHeadstagelNChipTempEnumNet Get↔
GroupChannelBitmaskHS1NChipTemp (
                                      uint32_t ChannelNumber,
                                       uint32_t virtualDevice )
11.8.2.23 GetGroupChannelBitmaskHS1Sidebands() [1/2] CMOSMeaHS1SidebandEnumNet GetGroup↔
ChannelBitmaskHS1Sidebands (
                                       uint32_t ChannelNumber )
\textbf{11.8.2.24} \quad \textbf{GetGroupChannelBitmaskHS1Sidebands() [2/2]} \quad \texttt{CMOSMeaHS1SidebandEnumNet} \quad \texttt{GetGroup} \leftarrow \texttt{CMOSMeaHS1SidebandEnumNet} \quad \texttt{CMOSMeaHS1Si
ChannelBitmaskHS1Sidebands (
                                       uint32_t ChannelNumber,
                                       uint32_t virtualDevice )
11.8.2.25 GetGroupChannelBitmaskHS1TriggerStatus() [1/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
GroupChannelBitmaskHS1TriggerStatus (
                                       uint32_t ChannelNumber )
11.8.2.26 GetGroupChannelBitmaskHS1TriggerStatus() [2/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
GroupChannelBitmaskHS1TriggerStatus (
                                       uint32_t ChannelNumber,
                                       uint32_t virtualDevice )
11.8.2.27 GetGroupChannelBitmaskIFDigChannels() [1/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
ChannelBitmaskIFDigChannels (
                                     uint32_t ChannelNumber )
11.8.2.28 GetGroupChannelBitmaskIFDigChannels() [2/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
ChannelBitmaskIFDigChannels (
                                       uint32_t ChannelNumber,
                                       uint32_t virtualDevice )
11.8.2.29 GetGroupChannelBitmaskInterfaceADC() [1/2] CMOSMeaInterfaceADCEnumNet GetGroup↔
ChannelBitmaskInterfaceADC (
                                       uint32_t ChannelNumber )
```

```
11.8.2.30 GetGroupChannelBitmaskInterfaceADC() [2/2] CMOSMeaInterfaceADCEnumNet GetGroup←
ChannelBitmaskInterfaceADC (
             uint32_t ChannelNumber,
             uint32_t virtualDevice )
11.8.2.31 GetGroupChannelBitmaskPacketFrameContext() [1/2] CMOSMeaPacketFrameContextGroup↔
EnumNet GetGroupChannelBitmaskPacketFrameContext (
              uint32_t ChannelNumber )
\textbf{11.8.2.32} \quad \textbf{GetGroupChannelBitmaskPacketFrameContext()} \  \  \texttt{[2/2]} \quad \texttt{CMOSMeaPacketFrameContextGroup} \leftarrow \\
{\tt EnumNet~GetGroupChannelBitmaskPacketFrameContext~(}
              uint32_t ChannelNumber,
              uint32_t virtualDevice )
11.8.2.33 GetGroupChannelBitmaskSTG1DACSignal() [1/2] CMOSMeaSTG1DACSignalEnumNet GetGroup←
ChannelBitmaskSTG1DACSignal (
              uint32_t ChannelNumber )
11.8.2.34 GetGroupChannelBitmaskSTG1DACSignal() [2/2] CMOSMeaSTG1DACSignalEnumNet GetGroup↔
ChannelBitmaskSTG1DACSignal (
             uint32_t ChannelNumber,
              uint32_t virtualDevice )
11.8.2.35 GetGroupDCOffset() [1/2] int32_t GetGroupDCOffset (
              DacqGroupChannelEnumNet GroupID )
\textbf{11.8.2.36} \quad \textbf{GetGroupDCOffset() [2/2]} \quad \texttt{int32\_t GetGroupDCOffset ()}
              DacqGroupChannelEnumNet GroupID,
              uint32_t virtualDevice )
11.8.2.37 GetGroupID() [1/2] DacqGroupChannelEnumNet GetGroupID (
              uint32_t Index )
```

```
11.8.2.38 GetGroupID() [2/2] DacqGroupChannelEnumNet GetGroupID (
             uint32_t Index,
             uint32_t virtualDevice )
11.8.2.39 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumNet GroupID )
11.8.2.40 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.8.2.41 GetGroupResolutionPerDigit() [1/2] uint32_t GetGroupResolutionPerDigit (
             DacqGroupChannelEnumNet GroupID )
11.8.2.42 GetGroupResolutionPerDigit() [2/2] uint32_t GetGroupResolutionPerDigit (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.8.2.43 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumNet GroupID )
11.8.2.44 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.8.2.45 GetGroupType() [1/2] DacqMeaGroupTypeEnumNet GetGroupType (
             {\tt DacqGroupChannelEnumNet}\ \textit{GroupID}\ )
11.8.2.46 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
```

```
11.8.2.47 GetGroupUnit() [1/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
             DacqGroupChannelEnumNet GroupID )
11.8.2.48 GetGroupUnit() [2/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.8.2.49 GetMaxNumOfColumns() uint32_t GetMaxNumOfColumns (
             uint32_t Samplerate )
\textbf{11.8.2.50} \quad \textbf{GetNeurochipMemoryData() [1/2]} \quad \texttt{uint} \ \texttt{32\_t GetNeurochipMemoryData} \ \ \textbf{(}
             uint16_t MemAddress )
11.8.2.51 GetNeurochipMemoryData() [2/2] array<uint32_t> ^ GetNeurochipMemoryData (
             uint16_t MemAddress,
             uint32_t ReqestLength )
11.8.2.52 GetNeurochipMemorySize() uint32_t GetNeurochipMemorySize ( )
11.8.2.53 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ( )
11.8.2.54 GetNumberOfSupportedGroups() [2/2] uint32_t GetNumberOfSupportedGroups (
             uint32_t virtualDevice )
11.8.2.55 GetSourceBulk() int32_t GetSourceBulk ( )
11.8.2.56 GetSourceDrain() int32_t GetSourceDrain ( )
```

```
11.8.2.57 GetSourceGate() int32_t GetSourceGate ( )
11.8.2.58 GetStimulusSites() List<int16_t> ^{\land} GetStimulusSites ( )
11.8.2.59 GetVDD3I() int32_t GetVDD3I ()
11.8.2.60 GetVDDI() int32_t GetVDDI ()
11.8.2.61 IsChipPowered() bool IsChipPowered ( )
\textbf{11.8.2.62} \quad \textbf{IsGateFloating()} \quad \texttt{bool IsGateFloating ()} \\
11.8.2.63 PowerChip() void PowerChip (
             bool on )
11.8.2.64 SetADCInputOffset() void SetADCInputOffset (
              int32_t offset )
11.8.2.65 SetBath() void SetBath (
              int32_t voltage )
11.8.2.66 SetBathMode() void SetBathMode (
              CMOSMeaBathModeEnumNet Mode )
```

```
11.8.2.67 SetGate() void SetGate (
             int32_t voltage )
11.8.2.68 SetGateFloating() void SetGateFloating ( )
11.8.2.69 SetGateToVOP() void SetGateToVOP ()
11.8.2.70 SetNeurochipMemoryData() [1/2] void SetNeurochipMemoryData (
             uint16_t MemAddress,
              array < uint32_t >^{\land} MemData)
11.8.2.71 SetNeurochipMemoryData() [2/2] void SetNeurochipMemoryData (
             uint16_t MemAddress,
             uint32_t MemData )
11.8.2.72 SetSourceBulk() void SetSourceBulk (
             int32_t voltage )
11.8.2.73 SetSourceDrain() void SetSourceDrain (
             int32_t voltage )
11.8.2.74 SetSourceGate() void SetSourceGate (
             int32_t voltage )
11.8.2.75 SetStimulusSites() void SetStimulusSites (
             List< int16_t >^{\land} SwitchPosition )
\textbf{11.8.2.76} \quad \textbf{UpdateTransistorVoltages()} \quad \texttt{void UpdateTransistorVoltages ()}
```


```
11.8.2.77 VOPSTimerSetResetTimes() [1/2] void VOPSTimerSetResetTimes ( uint32_t ResetTime, uint32_t IntervalTime )
```

11.8.2.78 VOPSTimerSetResetTimes() [2/2] void VOPSTimerSetResetTimes (uint32_t ResetTime,

uint32_t IntervalTime,
uint32_t HPFilterResetTime)

11.9 CCMOSMeaDeviceNet Class Reference

Inheritance diagram for CCMOSMeaDeviceNet:

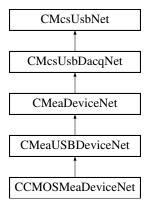

Classes

class CRegionOfInterestRect

Public Member Functions

- CCMOSMeaDeviceNet (void)
- ∼CCMOSMeaDeviceNet ()
- void SetBaseSamplerate (int BaseSamplerate)
- int GetBaseSamplerate ()
- array< int > ^ GetAvailableBaseSamplerates ()
- void SetRegionOfInterests (System::Collections::Generic::Dictionary< int, CRegionOfInterestRect^>^ rois)
- void UpdateChannelBlock (int queuesize, int threshold, int channels_in_block)
- System::Collections::Generic::Dictionary< int, array< array< int16_t >^> ^ GetCMOSDataDictionary (int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ GetChannelDataUI16 (DacqGroup← ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ GetChannelDatal16 (DacqGroup← ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ GetChannelDataUl32 (DacqGroup← ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames ret)
- System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ GetChannelDatal32 (DacqGroup← ChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Properties

```
• CCMOSMea_FunctionNet^ CMosMea [get]
```

```
• CStimulusFunctionNet^ Stimulus [get]
```

Additional Inherited Members

```
11.9.1 Constructor & Destructor Documentation
```

```
11.9.1.1 CCMOSMeaDeviceNet() CCMOSMeaDeviceNet (
            void )
11.9.1.2 ~CCMOSMeaDeviceNet() ~CCMOSMeaDeviceNet ()
11.9.2 Member Function Documentation
11.9.2.1 GetAvailableBaseSamplerates() array<int> ^ GetAvailableBaseSamplerates ()
11.9.2.2 GetBaseSamplerate() int GetBaseSamplerate ( )
11.9.2.3 GetChannelDatal16() System::Collections::Generic::Dictionary<int, array<int16_t>^> ^
GetChannelDataI16 (
            DacqGroupChannelEnumNet group,
            int frames,
            [System::Runtime::InteropServices::Out] int % frames_ret )
11.9.2.4 GetChannelDatal32() System::Collections::Generic::Dictionary<int, array<int32_t>^> ^
GetChannelDataI32 (
```

DacqGroupChannelEnumNet group,


[System::Runtime::InteropServices::Out] int % frames_ret)

int frames,

```
11.9.2.5 GetChannelDataUl16() System::Collections::Generic::Dictionary<int, array<uint16_t>^>
^ GetChannelDataUI16 (
             DacqGroupChannelEnumNet group,
             int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.9.2.6 GetChannelDataUl32() System::Collections::Generic::Dictionary<int, array<uint32_t>^>
^ GetChannelDataUI32 (
             DacqGroupChannelEnumNet group,
             int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.9.2.7 GetCMOSDataDictionary() System::Collections::Generic::Dictionary<int, array<array<int16↔
_{t}^{\wedge} \(^{\} \) \(^{\} \) GetCMOSDataDictionary (
             int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.9.2.8 SetBaseSamplerate() void SetBaseSamplerate (
             int BaseSamplerate )
11.9.2.9 SetRegionOfInterests() void SetRegionOfInterests (
             {\tt System::Collections::Generic::Dictionary<~int,~CRegionOfInterestRect^>^~\it rois~)}
11.9.2.10 UpdateChannelBlock() void UpdateChannelBlock (
             int queuesize,
             int threshold,
             int channels_in_block )
11.9.3 Property Documentation
11.9.3.1 CMosMea CCMOSMea_FunctionNet^ CMosMea [get]
11.9.3.2 Stimulus CStimulusFunctionNet^ Stimulus [get]
```

11.10 CCreateFilterNet Class Reference

Inheritance diagram for CCreateFilterNet:

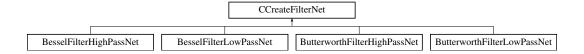

Public Member Functions

- · CCreateFilterNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)
- ∼CCreateFilterNet ()
- CFilterCoefficientsNet ^ GetBiQuad (int index)
- array< CFilterCoefficientsNet[^]> [^] GetBiQuads ()

Static Public Member Functions

- static int FindFilter (array< CFilterCoefficientsNet^>^ coef, array< CCreateFilterNet^>^ param)
- static int FindFilter (array< array< uint64_t >^> coef, array< CCreateFilterNet^> param, CFilterCoefficientsNet::s_FilterAttributesNet^ FiltAttr, bool DoMCSLegacyCompare)

Protected Member Functions

• CCreateFilterNet (int numCoefSets, CCreateFilter *pCreateFilter)

Properties

- int NumCoefSets [get]
- int Order [get]
- double SampleRate [get]
- double CutoffFrequency [get]
- double Scale [get]

11.10.1 Constructor & Destructor Documentation

11.10.1.1 CCreateFilterNet() [1/2] CCreateFilterNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

```
11.10.1.2 ~CCreateFilterNet() ~CCreateFilterNet ()
11.10.1.3 CCreateFilterNet() [2/2] CCreateFilterNet (
             int numCoefSets,
             CCreateFilter * pCreateFilter ) [protected]
11.10.2 Member Function Documentation
11.10.2.1 FindFilter() [1/2] static int FindFilter (
            array< uint64_t >^{\wedge}>^{\wedge} coef,
             array< CCreateFilterNet^>^ param,
            CFilterCoefficientsNet::s_FilterAttributesNet^ FiltAttr,
            bool DoMCSLegacyCompare ) [static]
11.10.2.2 FindFilter() [2/2] static int FindFilter (
            array< CFilterCoefficientsNet^>^ coef,
             array< CCreateFilterNet^>^ param ) [static]
11.10.2.3 GetBiQuad() CFilterCoefficientsNet ^ GetBiQuad (
            int index )
11.10.2.4 GetBiQuads() array<CFilterCoefficientsNet^> ^ GetBiQuads ()
11.10.3 Property Documentation
11.10.3.1 CutoffFrequency double CutoffFrequency [get]
11.10.3.2 NumCoefSets int NumCoefSets [get]
```


11.10.3.3 Order int Order [get]

11.10.3.4 SampleRate double SampleRate [get]

11.10.3.5 Scale double Scale [get]

11.11 CDacCalibrationFunctionNet Class Reference

Inheritance diagram for CDacCalibrationFunctionNet:

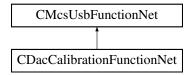

Public Member Functions

CDacCalibrationFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pDac
 — CalibrationFunctionPointerContainer)

Initializes a new instance of the CDacCalibrationFunctionNet class.

- CDacCalibrationFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CDacCalibrationFunctionNet ()
- !CDacCalibrationFunctionNet ()
- void SetDacOffset (uint16_t dacChannel, uint32_t offset)

Sets the offset of a DAC channel.

• uint32_t GetDacOffset (uint16_t dacChannel)

Gets the offset of a DAC channel.

void BurnDacOffset (uint16_t dacChannel)

Writes the offset of a DAC channel to permanent memory.

Additional Inherited Members

11.11.1 Detailed Description

11.11.2 Constructor & Destructor Documentation

```
11.11.2.1 CDacCalibrationFunctionNet() [1/2] CDacCalibrationFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pDacCalibrationFunctionPointerContainer)
```

Initializes a new instance of the CDacCalibrationFunctionNet class.

```
11.11.2.2 CDacCalibrationFunctionNet() [2/2] CDacCalibrationFunctionNet (
CMcsUsbNet^ mcsusb)
```

```
11.11.2.3 ~CDacCalibrationFunctionNet() virtual ~CDacCalibrationFunctionNet () [virtual]
```

```
11.11.2.4 "!CDacCalibrationFunctionNet() !CDacCalibrationFunctionNet ( )
```

11.11.3 Member Function Documentation

```
11.11.3.1 BurnDacOffset() void BurnDacOffset (
uint16_t dacChannel)
```

Writes the offset of a DAC channel to permanent memory.

Parameters

```
dacChannel The DAC channel number.
```

Gets the offset of a DAC channel.

Parameters

dacChannel	The DAC channel number.

Returns

The offset in digits.


Sets the offset of a DAC channel.

Parameters

dacChannel	The DAC channel number.
offset	The offset in digits.

11.12 CDacqGroupChannelGenericSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelGenericSelectionNet:

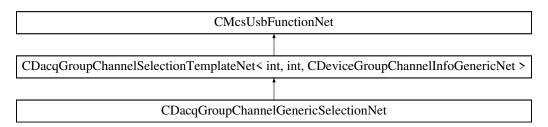

Public Member Functions

CDacqGroupChannelGenericSelectionNet (CMcsUsbNet[^] mcsusb)


Additional Inherited Members

11.12.1 Constructor & Destructor Documentation

```
11.12.1.1 CDacqGroupChannelGenericSelectionNet() CDacqGroupChannelGenericSelectionNet ( CMcsUsbNet^ mcsusb )
```

11.13 CDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelSelectionNet:

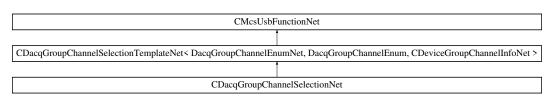

Public Member Functions

CDacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members

11.13.1 Constructor & Destructor Documentation

11.13.1.1 CDacqGroupChannelSelectionNet() CDacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)

11.14 CDacqGroupChannelSelectionTemplateNet < DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference

Inheritance diagram for CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet >:



Public Member Functions

- CDacqGroupChannelSelectionTemplateNet (CMcsUsbNet^{\(\Lambda\)} mcsusb)
- uint32_t GetNumberOfSupportedGroups ()
- uint32_t GetNumberOfSupportedGroups (uint32_t virtualDevice)
- DacqGroupChannelEnumTemplateNet GetGroupID (uint32_t Index)
- DacqGroupChannelEnumTemplateNet GetGroupID (uint32_t Index, uint32_t virtualDevice)
- uint32_t GetGroupNumberOfChannels (DacqGroupChannelEnumTemplateNet GroupID)
- DacqMeaGroupTypeEnumNet GetGroupType (DacqGroupChannelEnumTemplateNet GroupID)
- DacqMeaGroupTypeEnumNet GetGroupType (DacqGroupChannelEnumTemplateNet GroupID, uint32_

 t virtualDevice)
- void EnableChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, List< bool >[^] Enabled←
 ChannelsBitMap)
- void EnableChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, List< bool >[^] Enabled←
 ChannelsBitMap, uint32_t virtualDevice)
- $\bullet \ \ \, \mathsf{List} \! < \mathsf{bool} > \land \ \, \mathsf{GetEnabledChannelsInGroup} \, (\mathsf{DacqGroupChannelEnumTemplateNet} \, \mathsf{GroupID})$
- List< bool > [^] GetEnabledChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, uint32_← t virtualDevice)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumTemplateNet GroupID)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumTemplateNet GroupID, uint32_t virtual ← Device)
- List< CDeviceGroupChannelInfoTemplateNet^> ^ GetDeviceGroupChannelInfos ()
- $\bullet \ \, \mathsf{List} \! < \mathsf{CDeviceGroupChannelInfoTemplateNet}^{\wedge} \! > \ \, ^{\wedge} \ \, \mathsf{GetDeviceGroupChannelInfoS} \ \, (\mathsf{uint32_t} \ \mathsf{virtualDevice}) \\$

11.14.1 Constructor & Destructor Documentation

```
11.14.1.1 CDacqGroupChannelSelectionTemplateNet() CDacqGroupChannelSelectionTemplateNet (
CMcsUsbNet^ mcsusb )
```

11.14.2 Member Function Documentation

```
11.14.2.1 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (

DacqGroupChannelEnumTemplateNet GroupID,

List< bool >^ EnabledChannelsBitMap )
```

```
11.14.2.2 EnableChannelsInGroup() [2/2] void EnableChannelsInGroup (

DacqGroupChannelEnumTemplateNet GroupID,

List< bool >^ EnabledChannelsBitMap,

uint32_t virtualDevice )
```


```
11.14.2.5 GetEnabledChannelsInGroup() [1/2] List<bool> ^ GetEnabledChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID )
```

```
\textbf{11.14.2.7} \quad \textbf{GetGroupID()} \; \texttt{[1/2]} \quad \texttt{DacqGroupChannelEnumTemplateNet} \; \texttt{GetGroupID} \; \; \texttt{(}
               uint32_t Index )
11.14.2.8 GetGroupID() [2/2] DacqGroupChannelEnumTemplateNet GetGroupID (
               uint32_t Index,
               uint32_t virtualDevice )
\textbf{11.14.2.9} \quad \textbf{GetGroupNumberOfChannels() [1/2]} \quad \texttt{uint32\_t GetGroupNumberOfChannels ()}
               DacqGroupChannelEnumTemplateNet GroupID )
11.14.2.10 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
               DacqGroupChannelEnumTemplateNet GroupID,
               uint32_t virtualDevice )
11.14.2.11 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
               {\tt DacqGroupChannelEnumTemplateNet} \ \ {\tt GroupID} \ )
11.14.2.12 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
               DacqGroupChannelEnumTemplateNet GroupID,
               uint32_t virtualDevice )
\textbf{11.14.2.13} \quad \textbf{GetGroupType()} \; \texttt{[1/2]} \quad \texttt{DacqMeaGroupTypeEnumNet} \; \texttt{GetGroupType} \; \; \texttt{(}
               DacqGroupChannelEnumTemplateNet GroupID )
11.14.2.14 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
               DacqGroupChannelEnumTemplateNet GroupID,
               uint32_t virtualDevice )
11.14.2.15 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ( )
```

11.14.2.16 **GetNumberOfSupportedGroups()** [2/2] uint32_t GetNumberOfSupportedGroups (uint32_t virtualDevice)

11.15 CDeviceGroupChannelInfoGenericNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoGenericNet:

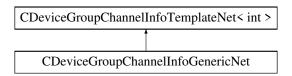

Public Member Functions


• CDeviceGroupChannelInfoGenericNet (int id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.15.1 Constructor & Destructor Documentation

11.16 CDeviceGroupChannelInfoNet Class Reference

 $Inheritance\ diagram\ for\ CDeviceGroupChannelInfoNet:$

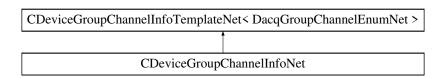

Public Member Functions

CDeviceGroupChannelInfoNet (DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

11.16.1 Constructor & Destructor Documentation

11.17 CDeviceGroupChannelInfoSCUNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoSCUNet:

```
CDeviceGroupChannelInfoTemplateNet < SCUDacqGroupChannelEnumNet >

CDeviceGroupChannelInfoSCUNet
```

Public Member Functions

Additional Inherited Members

11.17.1 Constructor & Destructor Documentation

11.18 CDeviceGroupChannelInfoTemplateNet < DacqGroupChannelEnumTemplateNet > Class Template Reference

Public Member Functions

• CDeviceGroupChannelInfoTemplateNet (DacqGroupChannelEnumTemplateNet id, int channels, DacqMea← GroupTypeEnumNet type)

Public Attributes

- DacqGroupChannelEnumTemplateNet GroupID
- int NumberOfChannels
- DacqMeaGroupTypeEnumNet GroupType

11.18.1 Constructor & Destructor Documentation

11.18.1.1 CDeviceGroupChannelInfoTemplateNet() CDeviceGroupChannelInfoTemplateNet (

```
DacqGroupChannelEnumTemplateNet id,
int channels,
DacqMeaGroupTypeEnumNet type )
```

11.18.2 Member Data Documentation


```
11.18.2.1 GroupID DacqGroupChannelEnumTemplateNet GroupID
```

11.18.2.2 GroupType DacqMeaGroupTypeEnumNet GroupType

11.18.2.3 NumberOfChannels int NumberOfChannels

11.19 CDeviceGroupChannelInfoW2100Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoW2100Net:

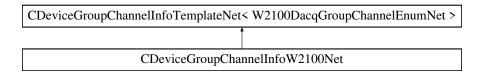

Public Member Functions


CDeviceGroupChannelInfoW2100Net (W2100DacqGroupChannelEnumNet id, int channels, DacqMea
 GroupTypeEnumNet type)

11.19.1 Constructor & Destructor Documentation

11.20 CDigOutStimulatorFunctionNet Class Reference

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

Inheritance diagram for CDigOutStimulatorFunctionNet:

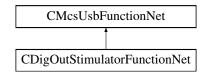

Public Member Functions

CDigOutStimulatorFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pDigOut
 — StimulatorFunctionPointerContainer)

Initializes a new instance of the CDigOutStimulatorFunctionNet class.

- CDigOutStimulatorFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CDigOutStimulatorFunctionNet ()
- !CDigOutStimulatorFunctionNet ()
- void ClearChannel (int32_t NrChannel)

clear stimulation pattern

- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareChannelData (array< int32_t >^ Amplitude, array< uint64 t >^ Duration)
- void SendChannelData (int32_t NrChannel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled)
- int32_t GetNumberOfChannels ()
- void SetGlobalRepeat (int32_t NrChannel, bool value)

set repeat whole stimulation pattern

bool GetGlobalRepeat (int32_t NrChannel)

get repeat whole stimulation pattern

• void SetStartTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)

sets start condition of digital out stimulator

• DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (int32_t NrChannel)

queries start condition of digital out stimulator

void SetStopTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)

sets stop condition of digital out stimulator

DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope (int32 t NrChannel)

queries stop condition of digital out stimulator

Additional Inherited Members

11.20.1 Detailed Description

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

11.20.2 Constructor & Destructor Documentation

```
11.20.2.1 CDigOutStimulatorFunctionNet() [1/2] CDigOutStimulatorFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pDigOutStimulatorFunctionPointerContainer )
```

Initializes a new instance of the CDigOutStimulatorFunctionNet class.

```
11.20.2.2 CDigOutStimulatorFunctionNet() [2/2] CDigOutStimulatorFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.20.2.3 ~CDigOutStimulatorFunctionNet() virtual ~CDigOutStimulatorFunctionNet () [virtual]

```
11.20.2.4 "!CDigOutStimulatorFunctionNet() !CDigOutStimulatorFunctionNet ( )
```

11.20.3 Member Function Documentation

```
11.20.3.1 ClearChannel() void ClearChannel (
    int32_t NrChannel )
```

clear stimulation pattern

Parameters

NrChannel

get repeat whole stimulation pattern

Parameters

NrChannel channel number

Returns

current value

11.20.3.3 GetNumberOfChannels() int32_t GetNumberOfChannels ()

Returns

11.20.3.4 GetStartTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (int32_t NrChannel)

queries start condition of digital out stimulator

Parameters

NrChannel channel number

Returns

start condition (rising or falling edge)

```
11.20.3.5 GetStopTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope ( int32_t NrChannel )
```

queries stop condition of digital out stimulator

Parameters

NrChannel channel number

Returns

stop condition (rising or falling edge)

Parameters

Amplitude	
Duration	

Returns

Parameters

NrChannel	
device_data_and_unrolled	

set repeat whole stimulation pattern

Parameters

NrChannel	channel number
value	new value

11.20.3.9 SetStartTriggerSlope() void SetStartTriggerSlope (

```
int32_t NrChannel,
DigitalStimulatorTriggerSlopeEnumNet Condition )
```

sets start condition of digital out stimulator

Parameters

NrChannel	channel number
Condition	start condition (rising or falling edge)

sets stop condition of digital out stimulator


Parameters

NrChannel	channel number
Condition	stop condition (rising or falling edge)

11.21 CEncapsulatorDeviceNet Class Reference

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

Inheritance diagram for CEncapsulatorDeviceNet:

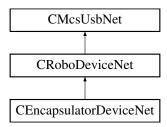

Public Member Functions

- CEncapsulatorDeviceNet (void)
- CRoboFluidDeviceNet ^ GetRoboFluidDevice ()

Additional Inherited Members

11.21.1 Detailed Description

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

11.21.2 Constructor & Destructor Documentation

```
11.21.2.1 CEncapsulatorDeviceNet() CEncapsulatorDeviceNet (
```


11.21.3 Member Function Documentation

```
11.21.3.1 GetRoboFluidDevice() CRoboFluidDeviceNet ^ GetRoboFluidDevice ( )
```

11.22 CExternDTesterDeviceNet Class Reference

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

Inheritance diagram for CExternDTesterDeviceNet:

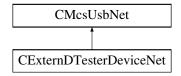

Public Member Functions

• CExternDTesterDeviceNet ()

Initializes a new instance of the CExternDTesterDeviceNet class.

- virtual \sim CExternDTesterDeviceNet ()
- !CExternDTesterDeviceNet ()
- array< uint8_t > $^{\land}$ Read (int configString_Length)

Reads the config string from the device.

String ^ Read2 ()

Reads the config string from the device.

void Write (array< uint8_t >^ configString)

Reads the config string from the device.

void Write2 (String[^] configString)

Reads the config string from the device.

Additional Inherited Members

11.22.1 Detailed Description

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

11.22.2 Constructor & Destructor Documentation

11.22.2.1 CExternDTesterDeviceNet() CExternDTesterDeviceNet ()

Initializes a new instance of the CExternDTesterDeviceNet class.

```
11.22.2.2 ~CExternDTesterDeviceNet() virtual ~CExternDTesterDeviceNet () [virtual]
```

```
11.22.2.3 "!CExternDTesterDeviceNet() !CExternDTesterDeviceNet ( )
```

11.22.3 Member Function Documentation

Reads the config string from the device.

Parameters

configString_Length	The maximal length of configString.
comigeting_congui	inc maxima length of comigoting.

Returns

The config string.

```
11.22.3.2 Read2() String ^{\land} Read2 ( )
```

Reads the config string from the device.

Returns

The config string.

Reads the config string from the device.

configString	The config string.
--------------	--------------------

```
11.22.3.4 Write2() void Write2 (
String^ configString)
```

Reads the config string from the device.

Parameters

configString	The config string.
--------------	--------------------

11.23 CFilterCoefficientsNet Class Reference

Classes

• struct s FilterAttributesNet

Public Member Functions

- CFilterCoefficientsNet ()
- CFilterCoefficientsNet (double b0, double b1, double b2, double a1, double a2)
- CFilterCoefficientsNet (double b0, double b1, double a1)
- CFilterCoefficientsNet (array< double >[∧] b, array< double >[∧] a)
- ∼CFilterCoefficientsNet ()
- bool IsEqual (array< uint64_t >^ coefficients, s_FilterAttributesNet^ FiltAttr)
- bool IsEqual (array< uint64_t >^ coefficients, s_FilterAttributesNet^ FiltAttr, bool DoMCSLegacyCompare)
- uint64_t GetUintB (int index, s_FilterAttributesNet^ FiltAttr)
- uint64_t GetUintA (int index, s_FilterAttributesNet[^] FiltAttr)

Properties

```
array< double > ^ A [get]array< double > ^ B [get]
```

11.23.1 Constructor & Destructor Documentation

11.23.1.1 CFilterCoefficientsNet() [1/4] CFilterCoefficientsNet ()

```
11.23.1.2 CFilterCoefficientsNet() [2/4] CFilterCoefficientsNet (
             double b0,
             double b1,
             double b2,
             double a1,
             double a2 )
11.23.1.3 CFilterCoefficientsNet() [3/4] CFilterCoefficientsNet (
             double b0,
             double b1,
             double a1 )
11.23.1.4 CFilterCoefficientsNet() [4/4] CFilterCoefficientsNet (
             array< double >^{\wedge} b,
             array< double >^{\wedge} a)
11.23.1.5 ~ CFilterCoefficientsNet() ~ CFilterCoefficientsNet ()
11.23.2 Member Function Documentation
11.23.2.1 GetUintA() uint64_t GetUintA (
             int index,
             s_FilterAttributesNet^ FiltAttr )
11.23.2.2 GetUintB() uint64_t GetUintB (
             int index,
             s_FilterAttributesNet^ FiltAttr )
11.23.2.3 IsEqual() [1/2] bool IsEqual (
             array< uint64_t >^{\wedge} coefficients,
             s_FilterAttributesNet^ FiltAttr )
```


11.23.3 Property Documentation

```
11.23.3.1 A array< double>^ A [get]
```

```
11.23.3.2 \mathbf{B} array< double>^ B [get]
```

11.24 CFilterConfigurationNet Class Reference

Inheritance diagram for CFilterConfigurationNet:

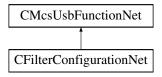

Public Member Functions

- CFilterConfigurationNet (CMcsUsbNet[^] mcsusb)
- void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet[^] Coefficients, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet[^] CoefficientsNet[^] CoefficientsNet[^] CoefficientsSet2, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32 t FilterNumber)
- void EraseFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)
- void SetHighpassFilterEnable (bool enable)
- bool GetHighpassFilterEnable ()
- void ResetHighpassFilter ()
- uint32_t GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, FilterAttribute

 EnumNet index)
- CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)


Additional Inherited Members

11.24.1 Constructor & Destructor Documentation

```
11.24.1.1 CFilterConfigurationNet() CFilterConfigurationNet (
             CMcsUsbNet^ mcsusb )
11.24.2 Member Function Documentation
11.24.2.1 EraseFilterParameterPermanent() void EraseFilterParameterPermanent (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber )
11.24.2.2 GetFilterAttributes() [1/2] CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilter←
Attributes (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber )
11.24.2.3 GetFilterAttributes() [2/2] uint32_t GetFilterAttributes (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber,
             FilterAttributeEnumNet index )
11.24.2.4 GetHighpassFilterEnable() bool GetHighpassFilterEnable ( )
11.24.2.5 ResetHighpassFilter() void ResetHighpassFilter ( )
11.24.2.6 SetFilterParameter() [1/2] void SetFilterParameter (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber,
             CFilterCoefficientsNet<sup>∧</sup> Coefficients,
             CFilterPropertyNet^ FilterProp )
11.24.2.7 SetFilterParameter() [2/2] void SetFilterParameter (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber,
             CFilterCoefficientsNet^ CoefficientsSet1,
             CFilterCoefficientsNet^ CoefficientsSet2,
             CFilterPropertyNet^ FilterProp )
```

11.25 CFilterConfigurationRegisterNet Class Reference

Inheritance diagram for CFilterConfigurationRegisterNet:

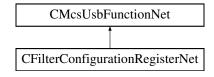

Public Member Functions

- CFilterConfigurationRegisterNet (CMcsUsbNet[^] mcsusb)
- void SetFilterParameter (uint32_t FilterCoefRegBase, CFilterCoefficientsNet[^] Coefficients, uint32_t Filter←
 InfoRegBase, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameter (uint32_t FilterCoefSet1RegBase, CFilterCoefficientsNet[^] CoefficientsSet1, uint32_t FilterCoefSet2RegBase, CFilterCoefficientsNet[^] CoefficientsSet2, uint32_t FilterInfoRegBase, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameterPermanent (uint32_t FilterCoefRegBase, uint32_t FilterCoefDmaReg, uint32_← t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBize)
- void SetFilterParameterPermanent (uint32_t FilterCoefSet1RegBase, uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2RegBase, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefDmaReg, uint32_t FilterInfoDmaReg, uint32_t E⇔
 EPROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)

Additional Inherited Members

11.25.1 Constructor & Destructor Documentation

```
11.25.1.1 CFilterConfigurationRegisterNet() CFilterConfigurationRegisterNet ( CMcsUsbNet^ mcsusb )
```

11.25.2 Member Function Documentation

```
11.25.2.1 EraseFilterParameterPermanent() [1/2] void EraseFilterParameterPermanent (
             uint32_t FilterCoefDmaReg,
             uint32_t FilterInfoDmaReg,
             uint32_t EEPROMBase,
             uint32_t \it EEPROMSize )
11.25.2.2 EraseFilterParameterPermanent() [2/2] void EraseFilterParameterPermanent (
             uint32_t FilterCoefSet1DmaReg,
             uint32_t FilterCoefSet2DmaReq,
             uint32_t FilterInfoDmaReg,
             uint32_t EEPROMBase,
             uint32_t EEPROMSize )
11.25.2.3 SetFilterParameter() [1/2] void SetFilterParameter (
             uint32_t FilterCoefRegBase,
             CFilterCoefficientsNet<sup>∧</sup> Coefficients,
             uint32_t FilterInfoRegBase,
             CFilterPropertyNet^ FilterProp )
11.25.2.4 SetFilterParameter() [2/2] void SetFilterParameter (
             uint32_t FilterCoefSet1RegBase,
             CFilterCoefficientsNet<sup>∧</sup> CoefficientsSet1,
             uint32_t FilterCoefSet2RegBase,
             CFilterCoefficientsNet^ CoefficientsSet2,
             uint32_t FilterInfoRegBase,
             CFilterPropertyNet^ FilterProp )
11.25.2.5 SetFilterParameterPermanent() [1/2] void SetFilterParameterPermanent (
             uint32_t FilterCoefRegBase,
             uint32_t FilterCoefDmaReg,
             uint32_t FilterInfoRegBase,
             uint32_t FilterInfoDmaReg,
             uint32_t EEPROMBase,
             uint32_t EEPROMSize )
11.25.2.6 SetFilterParameterPermanent() [2/2] void SetFilterParameterPermanent (
             uint32_t FilterCoefSet1RegBase,
             uint32_t FilterCoefSet1DmaReg,
             uint32_t FilterCoefSet2RegBase,
             uint32_t FilterCoefSet2DmaReg,
             uint32_t FilterInfoRegBase,
             uint32_t FilterInfoDmaReg,
             uint32_t EEPROMBase,
             uint32_t EEPROMSize )
```

11.26 CFilterPropertyNet Class Reference

Public Member Functions

- FamilyEnumNet FilterFamily, FilterTypeEnumNet FilterType, bool Active)
- ∼CFilterPropertyNet ()
- virtual System::String ^ ToString () override

Properties

```
• uint32_t CornerFrequencymHz [get]
• double CornerFrequency [get]
• uint32_t Order [get]
• FilterBandEnumNet FilterBand [get]
```

- FilterFamilyEnumNet FilterFamily [get]
- FilterTypeEnumNet FilterType [get]
- bool FilterActive [get]

11.26.1 Constructor & Destructor Documentation

```
11.26.1.1 CFilterPropertyNet() CFilterPropertyNet (
             uint32_t CornerFrequenzymHz,
             uint32_t Order,
             FilterBandEnumNet FilterBand,
             FilterFamilyEnumNet FilterFamily,
             FilterTypeEnumNet FilterType,
             bool Active )
```

```
11.26.1.2 ~CFilterPropertyNet() ~CFilterPropertyNet ( )
```

11.26.2 Member Function Documentation

```
11.26.2.1 ToString() virtual System::String ^ ToString ( ) [override], [virtual]
```

11.26.3 Property Documentation

11.26.3.1 CornerFrequency double CornerFrequency [get]

11.26.3.2 CornerFrequencymHz uint32_t CornerFrequencymHz [get]

11.26.3.3 FilterActive bool FilterActive [get]

11.26.3.4 FilterBand FilterBandEnumNet FilterBand [get]

11.26.3.5 FilterFamily FilterFamilyEnumNet FilterFamily [get]


11.26.3.6 FilterType FilterTypeEnumNet FilterType [get]

11.26.3.7 Order uint32_t Order [get]

11.27 CFluidControlDeviceNet Class Reference

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

Inheritance diagram for CFluidControlDeviceNet:

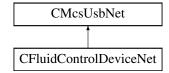

Public Member Functions

CFluidControlDeviceNet ()

Initialize a new instance of the CFluidControlDeviceNet class.

∼CFluidControlDeviceNet ()

Default destructor.

void SetValve (unsigned int value)

Open or Close valves.

· void SetSingleValve (unsigned short valve, unsigned short onoff)

Opens or Closes a valve.

void SetDigout (unsigned int value)

Define the pattern on the Digital Output.

void SetPWM (unsigned int channel, unsigned int value)

Sets the duty cycle of the PWM output.

void CalibrateThermocouple (unsigned int channel)

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

• void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)

Sets the proportinal constant for the Thermocouple.

unsigned int GetValve ()

Gets the state of the valves.

unsigned short GetSingleValve (unsigned short valve)

Gets the state of a valve.

unsigned int GetDigout ()

Gets the state of the digital output.

· unsigned int GetPWM (unsigned int channel)

Gets the state of the PWM output.

unsigned int GetAdc (unsigned int channel)

Reads an ADC Value.

• unsigned int GetDigin ()

Reads the digital input.

int GetThermocoupleTemperature (unsigned int channel)

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

• int GetReferenceTemperature (unsigned int channel)

Reads the reference temperature for the Thermocouple.

unsigned int GetThermocoupleCalibration (unsigned int channel)

Gets the calibration constant for the Thermocouple ADC.

unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)

Reads the proportional constant for the Thermocouple.

Properties

CMcsBus VoltageModeNet[^] McsBus VoltageMode [get]

Additional Inherited Members

11.27.1 Detailed Description

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

11.27.2 Constructor & Destructor Documentation

11.27.2.1 CFluidControlDeviceNet() CFluidControlDeviceNet ()

Initialize a new instance of the CFluidControlDeviceNet class.

11.27.2.2 ~CFluidControlDeviceNet() ~CFluidControlDeviceNet ()

Default destructor.

11.27.3 Member Function Documentation

11.27.3.1 CalibrateThermocouple() void CalibrateThermocouple (unsigned int channel)

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

Parameters

channel Thermocouple channel number.

11.27.3.2 **GetAdc()** unsigned int GetAdc (unsigned int *channel*)

Reads an ADC Value.

Parameters

channel	The ADC channel number to query.

Returns

The current ADC value.

11.27.3.3 **GetDigin()** unsigned int GetDigin ()

Reads the digital input.

Returns

The bit pattern of the state of the digital inputs.

11.27.3.4 GetDigout() unsigned int GetDigout ()

Gets the state of the digital output.

Returns

The current state of the digital outputs as a bit pattern.

11.27.3.5 **GetPWM()** unsigned int GetPWM (unsigned int *channel*)

Gets the state of the PWM output.

Returns

The current state of the PWM outputs duty cycle in permille.

11.27.3.6 GetReferenceTemperature() int GetReferenceTemperature (unsigned int *channel*)

Reads the reference temperature for the Thermocouple.

Parameters

Returns

The temperature from the Thermocouple in 1/100 °C.

11.27.3.7 GetSingleValve() unsigned short GetSingleValve (unsigned short *valve*)

Gets the state of a valve.

valve number of valve

Returns

state of the valve

11.27.3.8 **GetThermocoupleCalibration()** unsigned int GetThermocoupleCalibration (unsigned int *channel*)

Gets the calibration constant for the Thermocouple ADC.

Parameters

channel	Thermocouple channel number.
---------	------------------------------

Returns

The calibration constant for the Thermocouple ADC.

11.27.3.9 **GetThermocoupleNanovoltPerKelvin()** unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int *channel*)

Reads the proportional constant for the Thermocouple.

Parameters

channel	Thermocouple channel number.
---------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

11.27.3.10 **GetThermocoupleTemperature()** int GetThermocoupleTemperature (unsigned int *channel*)

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

channel	Thermocouple channel number.
---------	------------------------------

Returns

The temperature difference between both Thermocouple junctions in 1/100 $^\circ\! C.$

11.27.3.11 GetValve() unsigned int GetValve ()

Gets the state of the valves.

Returns

The current state of the valves as a bit pattern.

```
11.27.3.12 SetDigout() void SetDigout (
unsigned int value)
```

Define the pattern on the Digital Output.

Parameters

n the digital output.	value
-----------------------	-------

```
11.27.3.13 SetPWM() void SetPWM (
unsigned int channel,
unsigned int value)
```

Sets the duty cycle of the PWM output.

channel	PWM channel number.
value	duty cycle of the PWM output in permille.

```
11.27.3.14 SetSingleValve() void SetSingleValve (
    unsigned short valve,
    unsigned short onoff)
```

Opens or Closes a valve.

<i>valve</i> num	per of valve to be changed.
------------------	-----------------------------

Parameters

onoff open or close the valv

11.27.3.15 SetThermocoupleNanovoltPerKelvin() void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)

Sets the proportinal constant for the Thermocouple.

Parameters

channel	Thermocouple channel number.	
value	proportinal constant for the Thermocouple in Nanovolt per Kelvin.	

11.27.3.16 SetValve() void SetValve (unsigned int value)

Open or Close valves.

Parameters

value bit pattern of valves which should be open.


11.27.4 Property Documentation

11.27.4.1 McsBus_VoltageMode CMcsBus_VoltageModeNet^ McsBus_VoltageMode [get]

11.28 CFYIDeviceNet Class Reference

CFYIDeviceNet is the class to control the MCS FYI device

Inheritance diagram for CFYIDeviceNet:

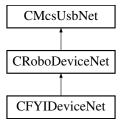

Public Member Functions

• CFYIDeviceNet (void)

Properties

- CRobo_FYITemp_FunctionNet^ FYITemp [get]
- CRobo_FYIProgram_FunctionNet^ FYIProgram [get]
- CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.28.1 Detailed Description

CFYIDeviceNet is the class to control the MCS FYI device

11.28.2 Constructor & Destructor Documentation

```
11.28.2.1 CFYIDeviceNet() CFYIDeviceNet ( void )
```

11.28.3 Property Documentation

11.28.3.1 FYIProgram CRobo_FYIProgram_FunctionNet^ FYIProgram [get]


11.28.3.2 FYITemp CRobo_FYITemp_FunctionNet^ FYITemp [get]

11.28.3.3 Sensor CMcsBus_SensorNet^ Sensor [get]

11.29 CGenericDevelopDeviceNet Class Reference

CGenericDevelopDeviceNet is the class to use during development of a new device.

Inheritance diagram for CGenericDevelopDeviceNet:

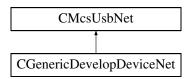

Public Member Functions

- CGenericDevelopDeviceNet (void)
 - Initialize a new instance of the CGenericDevelopDeviceNet class.
- ~CGenericDevelopDeviceNet (void)
- void SetValue (uint16_t value, uint16_t index)

Sets .

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index The index of the request.

- template<typename C >
 void SetBuffer (uint16_t value, uint16_t index, array< C >^ buffer)
- void SetUByteBuffer (uint16_t value, uint16_t index, array< unsigned char >^ buffer)
 Sends an array of type unsigned char to the device.

Parameters

value	The value of the request.

index	The index of the request.
-------	---------------------------

buffer	The buffer to send.
--------	---------------------

• void SetByteBuffer (uint16_t value, uint16_t Index, array< char $>^{\land}$ buffer) Sends an array of type char to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

Index	The index of the request.
-------	---------------------------

Parameters

• void SetUShortBuffer (uint16_t value, uint16_t index, array< unsigned short >^ buffer)

Sends an array of type unsigned short to the device.

Parameters

value	The value of the request.

Parameters

index	The index of the request.

Parameters

• void SetShortBuffer (uint16_t value, uint16_t index, array< short $>^{\land}$ buffer) Sends an array of type short to the device.

value	The value of the request.

index The index of the request.

Parameters

buffer	The buffer to send.

• void SetUIntBuffer (uint16_t value, uint16_t index, array< unsigned int >^ buffer)

Sends an array of unsigned int to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

buffer The buffer to send.

• void SetIntBuffer (uint16_t value, uint16_t index, array< int >^ buffer)

Sends an array of type int to the device.

Parameters

value	The value of the request.

Parameters

index	The index of the request.
-------	---------------------------

buffer The buffer to ser	d.
--------------------------	----

template<typename C >

array< C > $^{\land}$ GetBuffer (uint16_t value, uint16_t index, int size)

array< unsigned char > ^ GetUByteBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned char from the device.

Parameters

value The va	lue of the request.
--------------	---------------------

Parameters

index The index of the request.	
---------------------------------	--

Parameters

Returns

The array of data from the device.

array< char > ^ GetByteBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type char from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index Th	ne index of the request.
----------	--------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

array< unsigned short > ^ GetUShortBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned short from the device.

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

array< short > ^ GetShortBuffer (uint16_t value, uint16_t index, int size)
 Gets an array of type short from the device.

Parameters

value The value of the	e request.
------------------------	------------

Parameters

index The index of the re

Parameters

Returns

The array of data from the device.

• array< unsigned int > ^ GetUIntBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned int from the device.

value	The value of the request.
-------	---------------------------

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

- array< int > ^ GetIntBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type int from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size The	size of the array.
----------	--------------------

Returns

The array of data from the device.

• IntPtr OpenPipe (uint8_t endpointAddress)

Open a Pipe to an USB Endpoint.

Parameters

endpointAddress	The Endpoint Number.
-----------------	----------------------

Returns

A handle to the endpoint.

• void ClosePipe (IntPtr pipeHandle)

Close a Pipe to an USB Endpoint.

pipeHandle	The endpoint handle.
------------	----------------------

void ResetPipe (IntPtr pipeHandle)

Reset a Pipe to an USB Endpoint.

Parameters

pipeHandle The endpoint handle.

template < typename C >
 array < C > ^ ReadPipe (IntPtr pipeHandle, uint32_t size)
 Read data from an USB Endpoint.

Parameters

pipeHandle The endpoint handle.

Parameters

size Number of elements to read.

Returns

An array of data read.

• template<typename C > void WritePipe (IntPtr pipeHandle, array< C > $^{\wedge}$ buffer)

Write data to an USB Endpoint.

Parameters

pipeHandle The endpoint handle.

Parameters

buffer An array of data to write.

Additional Inherited Members

11.29.1 Detailed Description

CGenericDevelopDeviceNet is the class to use during development of a new device.

11.29.2 Constructor & Destructor Documentation

```
11.29.2.1 CGenericDevelopDeviceNet() CGenericDevelopDeviceNet (
```

Initialize a new instance of the CGenericDevelopDeviceNet class.

```
11.29.2.2 ~CGenericDevelopDeviceNet() ~CGenericDevelopDeviceNet ( void )
```

11.29.3 Member Function Documentation

Close a Pipe to an USB Endpoint.

Parameters

```
pipeHandle The endpoint handle.
```

Gets an array of type char from the device.

<i>value</i> T	he value of the request.
----------------	--------------------------

Parameters

index The index of the reque	st.
------------------------------	-----

Parameters

size The size of	of the array.
------------------	---------------

Returns

The array of data from the device.

Gets an array of type int from the device.

Parameters

value	The value of the request.
value	The value of the request.

index	The index of the request.

Returns

The array of data from the device.

Gets an array of type short from the device.

Parameters

	value	The value of the request.
--	-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

Gets an array of type unsigned char from the device.

<i>value</i> T	he value of the request.
----------------	--------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size The size of	of the array.
------------------	---------------

Returns

The array of data from the device.

Gets an array of type unsigned int from the device.

Parameters

index	The index of the request.

size The size of the arra	ay.
---------------------------	-----

Returns

The array of data from the device.

Gets an array of type unsigned short from the device.

Parameters

value The value of the reque	st.
------------------------------	-----

Parameters

index	The index of the request.
-------	---------------------------

Parameters

```
size The size of the array.
```

Returns

The array of data from the device.

```
11.29.3.9 OpenPipe() IntPtr OpenPipe ( uint8_t endpointAddress )
```

Open a Pipe to an USB Endpoint.

Returns

A handle to the endpoint.

Read data from an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.
------------	----------------------

Parameters

Returns

An array of data read.

Reset a Pipe to an USB Endpoint.

pipeHandle	The endpoint handle.

Sends an array of type char to the device.

Parameters

value The value of the requ	est.
-----------------------------	------

Parameters

Parameters

```
buffer The buffer to send.
```

Sends an array of type int to the device.

value	The value of the request.
-------	---------------------------

index	The index of the request.
-------	---------------------------

Parameters

buffer The buffer to send.

uint16_t *index*, array< short $>^{\land}$ buffer)

Sends an array of type short to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

Parameters

buffer The buffer to send.

11.29.3.16 SetUByteBuffer() void SetUByteBuffer (uint16_t value,

```
\label{eq:continuous} \mbox{uint16\_t } \mbox{index,} \\ \mbox{array< unsigned char } >^{\wedge} \mbox{buffer )}
```

Sends an array of type unsigned char to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

```
buffer The buffer to send.
```


Sends an array of unsigned int to the device.

Parameters

value The value of the request.	
---------------------------------	--

index The index of the request.	
---------------------------------	--

buffer The buffer to send.

Sends an array of type unsigned short to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

buffer	The buffer to send.
--------	---------------------

Sets.

Parameters

value The value of the request.

index The index of the i	request.
--------------------------	----------

Write data to an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.
------------	----------------------


Parameters

buffer	An array of data to write.
--------	----------------------------

11.30 CGilsonDeviceNet Class Reference

CGilsonDeviceNet is the class to control a Gilson device.

Inheritance diagram for CGilsonDeviceNet:

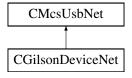

Public Member Functions

- CGilsonDeviceNet (void)
 - Initialize a new instance of the CGilsonDeviceNet class.
- ∼CGilsonDeviceNet (void)
- void ConnectSlave (byte ID)
- void SendImmediate (wchar_t command)
- String \(^\) SendImmediateGetResponse (wchar_t command)
- void SendBuffered (String[^] command)
- String ^ GetLastAnswer ()

Protected Attributes

CGilsonDevice * m_pGilsonDevice

Additional Inherited Members

11.30.1 Detailed Description

CGilsonDeviceNet is the class to control a Gilson device.

11.30.2 Constructor & Destructor Documentation

```
11.30.2.1 CGilsonDeviceNet() CGilsonDeviceNet (
void )
```

Initialize a new instance of the CGilsonDeviceNet class.

```
11.30.2.2 ~CGilsonDeviceNet() ~CGilsonDeviceNet (
```

11.30.3 Member Function Documentation

```
11.30.3.1 ConnectSlave() void ConnectSlave ( byte ID )
```

```
11.30.3.2 GetLastAnswer() String ^ GetLastAnswer ( )
```

```
11.30.3.3 SendBuffered() void SendBuffered ( String^{\wedge} command )
```

11.30.3.5 SendImmediateGetResponse() String ^ SendImmediateGetResponse (wchar_t command)


11.30.4 Member Data Documentation

11.30.4.1 m_pGilsonDevice CGilsonDevice* m_pGilsonDevice [protected]

11.31 CGrapheneFunctionNet Class Reference

CGrapheneFunctionNet is the class to control the TEER device

Inheritance diagram for CGrapheneFunctionNet:

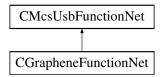

Public Member Functions

CGrapheneFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pGraphene ← FunctionPointerContainer)

Initializes a new instance of the CGrapheneFunctionNet class.

- CGrapheneFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CGrapheneFunctionNet ()
- !CGrapheneFunctionNet ()
- void GetVdVsDAC ([System::Runtime::InteropServices::Out]int16_t% Vd, [System::Runtime::Interop←
 Services::Out]int16_t% Vs)
- void GetVdVsDAC (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% Vd, [System::
 Runtime::InteropServices::Out]int16_t% Vs)
- void SetVdVsDAC (int16_t Vd, int16_t Vs)
- void SetVdVsDAC (uint32_t Headstage, int16_t Vd, int16_t Vs)
- bool GetVoltageReached ()
- bool GetVoltageReached (uint32_t Headstage)
- int32_t GetVoltageRange ()
- int32_t GetVoltageRange (uint32_t Headstage)
- void SetVoltageRange (int32 t range)
- void SetVoltageRange (uint32 t Headstage, int32 t range)
- int32_t GetVoltageResolution ()
- int32_t GetVoltageResolution (uint32_t Headstage)
- void SetVoltageResolution (int32_t resolution)
- void SetVoltageResolution (uint32_t Headstage, int32_t resolution)
- void GetDACOffset ([System::Runtime::InteropServices::Out]int16_t% offset_vd, [System::Runtime::
 — InteropServices::Out]int16 t% offset vs)
- void GetDACOffset (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% offset_vd, [System::Runtime::InteropServices::Out]int16_t% offset_vs)
- void SetDACOffset (int16_t offset_vd, int16_t offset_vs)

- void SetDACOffset (uint32_t Headstage, int16_t offset_vd, int16_t offset_vs)
- void GetVdVs ([System::Runtime::InteropServices::Out]int32_t% Vd, [System::Runtime::InteropServices::
 — Out]int32_t% Vs)
- void GetVdVs (uint32_t Headstage, [System::Runtime::InteropServices::Out]int32_t% Vd, [System::←
 Runtime::InteropServices::Out]int32_t% Vs)
- void SetVdVs (int32_t Vd, int32_t Vs)
- void SetVdVs (uint32 t Headstage, int32 t Vd, int32 t Vs)
- void SetVgs (int32_t Vgs)
- void SetVgs (uint32_t Headstage, int32_t Vgs)
- void SetVds (int32_t Vds)
- void SetVds (uint32_t Headstage, int32_t Vds)

Additional Inherited Members

11.31.1 Detailed Description

CGrapheneFunctionNet is the class to control the TEER device

11.31.2 Constructor & Destructor Documentation

```
11.31.2.1 CGrapheneFunctionNet() [1/2] CGrapheneFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pGrapheneFunctionPointerContainer)
```

Initializes a new instance of the CGrapheneFunctionNet class.

```
11.31.2.2 CGrapheneFunctionNet() [2/2] CGrapheneFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
\textbf{11.31.2.3} \quad \sim \textbf{CGrapheneFunctionNet()} \quad \text{virtual} \quad \sim \textbf{CGrapheneFunctionNet ()} \quad [\text{virtual}]
```

```
11.31.2.4 "!CGrapheneFunctionNet() !CGrapheneFunctionNet ( )
```

11.31.3 Member Function Documentation

offset_vd	
offset_vs	

11.31.3.2 GetDACOffset() [2/2] void GetDACOffset (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int16_t% offset_vd,
[System::Runtime::InteropServices::Out] int16_t% offset_vs )
```

Parameters

Headstage	The headstage to query.
offset_vd	
offset_vs	

11.31.3.3 GetVdVs() [1/2] void GetVdVs (

```
[System::Runtime::InteropServices::Out] int32_t% Vd, [System::Runtime::InteropServices::Out] int32_t% Vs)
```

Parameters

Vd	
Vs	

11.31.3.4 GetVdVs() [2/2] void GetVdVs (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int32_t% Vd,
[System::Runtime::InteropServices::Out] int32_t% Vs )
```

Parameters

Headstage	The headstage to query.
Vd	
Vs	

11.31.3.5 GetVdVsDAC() [1/2] void GetVdVsDAC (

```
[System::Runtime::InteropServices::Out] int16_t% Vd, [System::Runtime::InteropServices::Out] int16_t% Vs)
```

Vd	
Vs	

11.31.3.6 GetVdVsDAC() [2/2] void GetVdVsDAC (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int16_t% Vd,
[System::Runtime::InteropServices::Out] int16_t% Vs )
```

Parameters

Headstage	The headstage to query.
Vd	
Vs	

11.31.3.7 GetVoltageRange() [1/2] int32_t GetVoltageRange ()

Returns

11.31.3.8 GetVoltageRange() [2/2] int32_t GetVoltageRange (uint32_t Headstage)

Parameters

Headstage	The headstage to query.

Returns

11.31.3.9 GetVoltageReached() [1/2] bool GetVoltageReached ()

Returns

```
11.31.3.10 GetVoltageReached() [2/2] bool GetVoltageReached ( uint32_t Headstage )
```

Headstage The headstage to query.

Returns

11.31.3.11 GetVoltageResolution() [1/2] int32_t GetVoltageResolution ()

Returns

```
11.31.3.12 GetVoltageResolution() [2/2] int32_t GetVoltageResolution ( uint32_t Headstage )
```

Parameters

Headstage The headstage to query.

Returns

Parameters

offset_vd offset_vs

```
11.31.3.14 SetDACOffset() [2/2] void SetDACOffset (
```

uint32_t Headstage,

```
int16_t offset_vd,
int16_t offset_vs )
```

Headstage	The headstage to query.
offset_vd	
offset_vs	

```
11.31.3.15 SetVds() [1/2] void SetVds ( int32_t Vds )
```

Parameters

Vds

Parameters

Headstage	The headstage to query.
Vds	

Vd	
Vs	

Headstage	The headstage to query.
Vd	
Vs	

Parameters

Vd	
Vs	

Parameters

Headstage	The headstage to query.
Vd	
Vs	

```
11.31.3.21 SetVgs() [1/2] void SetVgs ( int32_t Vgs )
```

Parameters

Vgs

Headstage	The headstage to query.
Vgs	

Parameters

range

Parameters

Headstage	The headstage to query.
range	

Parameters

resolution

11.31.3.26 SetVoltageResolution() [2/2] void SetVoltageResolution (uint32_t Headstage, int32_t resolution)


Parameters

Headstage	The headstage to query.
resolution	

11.32 CHiClampDeviceNet Class Reference

CHiClampDeviceNet is the to control the MCS HiClamp device

Inheritance diagram for CHiClampDeviceNet:

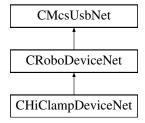

Public Member Functions

• CHiClampDeviceNet (void)

Properties

CRoboDacqNet[^] RoboDacq [get]

Additional Inherited Members

11.32.1 Detailed Description

CHiClampDeviceNet is the to control the MCS HiClamp device

11.32.2 Constructor & Destructor Documentation


```
11.32.2.1 CHiClampDeviceNet() CHiClampDeviceNet ( void )
```

11.32.3 Property Documentation

11.32.3.1 RoboDacq CRoboDacqNet^ RoboDacq [get]

11.33 CHLADacqNet Class Reference

Inheritance diagram for CHLADacqNet:

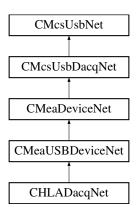

Public Member Functions

CHLADacqNet (void)

Additional Inherited Members


11.33.1 Constructor & Destructor Documentation

```
11.33.1.1 CHLADacqNet() CHLADacqNet (
void )
```

11.34 CHLADeviceNet Class Reference

CHLADeviceNet is the to control the MCS HLA device

Inheritance diagram for CHLADeviceNet:

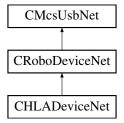

Public Member Functions

• CHLADeviceNet (void)

Properties

- CHLADacqNet^ HLADacq [get]
- CSerialPortNet^ SerialPort [get]

Additional Inherited Members

11.34.1 Detailed Description

CHLADeviceNet is the to control the MCS HLA device

11.34.2 Constructor & Destructor Documentation

```
11.34.2.1 CHLADeviceNet() CHLADeviceNet (
```

11.34.3 Property Documentation

```
11.34.3.1 HLADacq CHLADacqNet^ HLADacq [get]
```

```
11.34.3.2 SerialPort CSerialPortNet^ SerialPort [get]
```

11.35 CMcsUsbDacqNet::CHWInfo Class Reference

Class to provide hardware information about the device.

Classes

· class CVoltageRangeInfoNet

Public Member Functions

- CHWInfo (CMcsUsbDacqNet[^] device)

Get the number of analog channels the device supports.

Get the number of digital channels the device supports.

· virtual bool IsDigitalChannelDedicated ()

Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

- virtual uint32_t GetAvailableSampleRates ([System::Runtime::InteropServices::Out]System::Collections::
 Generic::List< int32_t >^% sampleRates)
- virtual uint32_t GetAvailableVoltageRangesInMicroVolt ([System::Runtime::InteropServices::Out]System::
 Collections::Generic::List< int32_t >^% voltageRanges)
- virtual uint32_t GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt ([System::Runtime::Interop← Services::Out]System::Collections::Generic::List< CVoltageRangeInfoNet^>^% voltageRanges)

11.35.1 Detailed Description

Class to provide hardware information about the device.

11.35.2 Constructor & Destructor Documentation

```
11.35.2.1 CHWInfo() CHWInfo (

CMcsUsbDacqNet^ device)
```

11.35.3 Member Function Documentation

```
11.35.3.4 GetNumberOfHWADCChannels() virtual uint32_t GetNumberOfHWADCChannels (
[System::Runtime::InteropServices::Out] int% numberOfChannels ) [virtual]
```

Get the number of analog channels the device supports.

Parameters

hberOfChannels Number of analog channels the device supports.

Returns

Error Status. 0 on success.

11.35.3.5 GetNumberOfHWDigitalChannels() virtual uint32_t GetNumberOfHWDigitalChannels ([System::Runtime::InteropServices::Out] int% numberOfChannels) [virtual]

Get the number of digital channels the device supports.

Parameters

numberOfChannels	Number of digital channels the device supports.
------------------	---

Returns

Error Status. 0 on success.

11.35.3.6 IsDigitalChannelDedicated() virtual bool IsDigitalChannelDedicated () [virtual]


Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

Returns

false when the digital channel replaces an analog channel when enabled, true when the digital channels is appended to the analog channels when enabled.

11.36 CIntanMea_FunctionNet Class Reference

Inheritance diagram for CIntanMea_FunctionNet:

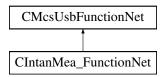

Public Member Functions

- CIntanMea_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] intalMea_Function
 —
 PointerContainer)
- CIntanMea FunctionNet (CMcsUsbNet[∧] mcsusb)
- int GetUpperFrequencyByIndex (unsigned short index)
- int GetLowerFrequencyByIndex (unsigned short index)
- int64_t GetDSPHighPassByIndex (unsigned short index)
- int GetIntanRegister (unsigned short chip, unsigned short registernumber)
- int GetImpedanceResult (unsigned short channel)
- void SetBandwidthByIndex (int upper_index, int lower_index)
- void SetDSPHighPassByIndex (int index)
- void AmplifierSettle ()
- void SetIntanRegister (unsigned short register_number, int value)
- void SetDiagnosticMode (unsigned char onoff)
- void BeginImpedanceCheck (array< int >^ config_values)

Additional Inherited Members

11.36.1 Constructor & Destructor Documentation


```
11.36.1.1 CIntanMea_FunctionNet() [1/2] CIntanMea_FunctionNet (
             CMcsUsbNet^ mcsusb,
             CMcsUsbFunctionPointerContainer^ intalMea_FunctionPointerContainer )
11.36.1.2 CIntanMea_FunctionNet() [2/2] CIntanMea_FunctionNet (
             CMcsUsbNet^ mcsusb )
11.36.2 Member Function Documentation
11.36.2.1 AmplifierSettle() void AmplifierSettle ( )
11.36.2.2 BeginImpedanceCheck() void BeginImpedanceCheck (
             array < int >^{\land} config\_values)
11.36.2.3 GetDSPHighPassByIndex() int64_t GetDSPHighPassByIndex (
             unsigned short index)
11.36.2.4 GetImpedanceResult() int GetImpedanceResult (
             unsigned short channel )
11.36.2.5 GetIntanRegister() int GetIntanRegister (
             unsigned short chip,
             unsigned short registernumber )
```

```
\textbf{11.36.2.6} \quad \textbf{GetLowerFrequencyByIndex()} \quad \texttt{int GetLowerFrequencyByIndex} \ \ \textbf{(}
              unsigned short index)
11.36.2.7 GetUpperFrequencyByIndex() int GetUpperFrequencyByIndex (
              unsigned short index)
11.36.2.8 SetBandwidthByIndex() void SetBandwidthByIndex (
              int upper_index,
              int lower_index )
11.36.2.9 SetDiagnosticMode() void SetDiagnosticMode (
              unsigned char onoff )
11.36.2.10 SetDSPHighPassByIndex() void SetDSPHighPassByIndex (
              int index )
11.36.2.11 SetIntanRegister() void SetIntanRegister (
              unsigned short register_number,
              int value )
```

11.37 CInterfaceboard2FunctionNet Class Reference

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboard2FunctionNet:

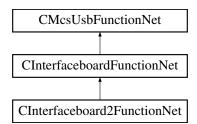

Public Member Functions

CInterfaceboard2FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pInterfaceboard2←
 FunctionPointerContainer)

Initializes a new instance of the CInterfaceboard2FunctionNet class.

- CInterfaceboard2FunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CInterfaceboard2FunctionNet ()
- !CInterfaceboard2FunctionNet ()
- void SetloVoltage (IoVoltageEnumNet ioVoltage)

Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

IoVoltageEnumNet GetIoVoltage ()

Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

Additional Inherited Members

11.37.1 Detailed Description

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

11.37.2 Constructor & Destructor Documentation

```
11.37.2.1 CInterfaceboard2FunctionNet() [1/2] CInterfaceboard2FunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pInterfaceboard2FunctionPointerContainer)
```

Initializes a new instance of the CInterfaceboard2FunctionNet class.

```
11.37.2.2 CInterfaceboard2FunctionNet() [2/2] CInterfaceboard2FunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.37.2.3 ~CInterfaceboard2FunctionNet() virtual ~CInterfaceboard2FunctionNet ( ) [virtual]
```

11.37.2.4 "!CInterfaceboard2FunctionNet() !CInterfaceboard2FunctionNet ()

11.37.3 Member Function Documentation

11.37.3.1 GetloVoltage() IoVoltageEnumNet GetIoVoltage ()

Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

Returns

Enum for the IO Voltage (3.3V or 5.0V).

Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.


Parameters

ioVoltage Enum for the I/O Voltage (3.3V or 5.0V).

11.38 CInterfaceboardFunctionNet Class Reference

CInterfaceboardFunctionNet is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboardFunctionNet:

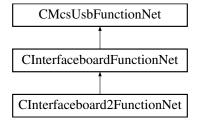

Public Member Functions

CInterfaceboardFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pInterfaceboard
 —
 FunctionPointerContainer)

Initializes a new instance of the CInterfaceboardFunctionNet class.

- CInterfaceboardFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CInterfaceboardFunctionNet ()
- !CInterfaceboardFunctionNet ()
- void SetCardinalDacqSamplerate (uint32_t samplerate)

Sets the fundamental/cardinal data aguisition samplerate of the Interfaceboard, default is 50 kHz

• uint32_t GetCardinalDacqSamplerate ()

Gets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz

void SetCardinalStgOutputrate (uint32_t outputrate)

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

uint32 t GetCardinalStgOutputrate ()

Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Additional Inherited Members

11.38.1 Detailed Description

CInterfaceboardFunctionNet is the class to control the Interfaceboard

11.38.2 Constructor & Destructor Documentation

Initializes a new instance of the CInterfaceboardFunctionNet class.

```
11.38.2.2 CInterfaceboardFunctionNet() [2/2] CInterfaceboardFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.38.2.3 ~CInterfaceboardFunctionNet() virtual ~CInterfaceboardFunctionNet () [virtual]
```

```
11.38.2.4 "!CInterfaceboardFunctionNet() !CInterfaceboardFunctionNet ()
```

11.38.3 Member Function Documentation

```
11.38.3.1 GetCardinalDacqSamplerate() uint32_t GetCardinalDacqSamplerate ( )
```

Gets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz Returns

The samplerate in Hz.

$\textbf{11.38.3.2} \quad \textbf{GetCardinalStgOutputrate()} \quad \texttt{uint32_t GetCardinalStgOutputrate ()} \\$

Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz $\,$

Returns

The output rate in Hz.

```
11.38.3.3 SetCardinalDacqSamplerate() void SetCardinalDacqSamplerate ( uint32_t samplerate )
```

Sets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz

samplerate

11.38.3.4 SetCardinalStgOutputrate() void SetCardinalStgOutputrate (uint32_t outputrate)

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz


Parameters

outputrate	The output rate in Hz.
------------	------------------------

11.39 CLIH3DeviceNet Class Reference

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

Inheritance diagram for CLIH3DeviceNet:

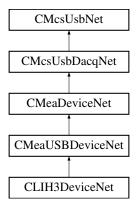

Public Member Functions

· CLIH3DeviceNet ()

Initializes a new instance of the CLIH3DeviceNet class.

- virtual ∼CLIH3DeviceNet ()
- !CLIH3DeviceNet ()
- void DummyCommand (uint32_t dummyParameter)

Dummy command to show how to use the DLL.

void SetEEpromPage (uint32_t EEpromStartAddress, array < int8_t >^ EEpromData, LIH30_EPC10_Bus
 _EnumNet epc10bus)

Writes into EEprom on the EPC10 EEPROM

array< int8_t > ^ GetEEpromPage (uint32_t EEpromStartAddress, int EEpromData_Length, LIH30_EP
 — C10_Bus_EnumNet epc10bus)

Reads the requested amount of EEprom byte from the EPC10 EEPROM

void SetSampleInterval (uint32_t SampleInterval)

Sets the Sample Interval for the DACQ and Stimulation

uint32 t GetSampleInterval ()

Gets the Sample Interval for the DACQ and Stimulation

void SetAdcOffset (LIH30 ADC Channel EnumNet AdcChannel, int32 t Offset)

Sets the ADC offset of the DACQ for a single channel

int32_t GetAdcOffset (LIH30_ADC_Channel_EnumNet AdcChannel)

Gets the ADC offset of the DACQ for a single channel

void SetAdcOffsetPermanent (LIH30_ADC_Channel_EnumNet AdcChannel)

Writes the ADC offset of the DACQ for a single channel to permanent EEProm memory

void ErasePermanentAdcOffset (LIH30_ADC_Channel_EnumNet AdcChannel)

Delets the ADC offset of the DACQ for a single channel in permanent EEProm memory

uint32_t ReadClipping (LIH30_EPC10_Bus_EnumNet epc10bus)

Gets the clipping information

void SetDigOutState (uint16_t DigOutState)

Writes to the LIH30 digital output

uint16_t GetDigInState ()

Reads from the LIH30 digital input

• void SendCommand (LIH30_EPC10_Bus_EnumNet epc10bus, uint16_t Command)

Send command to the EPC10

uint16_t GetDacqRunStatus ()

Gets the data acquisition running status

void SetDacUseIdleValue (uint32 t DacChannel, bool UseIdle)

Sets if the DAC Idle value is used after stimulation

bool GetDacUseIdleValue (uint32_t DacChannel)

Gets if the DAC Idle value is used after stimulation

void SetDacIdleValue (uint32_t DacChannel, int32_t IdleValue)

Sets the DAC Idle value

int32_t GetDacIdleValue (uint32_t DacChannel)

Gets the DAC Idle value

void EnableUserTrigger (bool enable)

Enables the User Trigger

bool IsUserTriggerEnabled ()

Is the User Trigger enabled

void SetDacOffset (LIH30 DAC Channel EnumNet DacChannel, int32 t Offset)

Sets the offset of a DAC channel.

int32_t GetDacOffset (LIH30_DAC_Channel_EnumNet DacChannel)

Gets the offset of a DAC channel.

void SetDacOffsetPermanent (LIH30_DAC_Channel_EnumNet DacChannel)

Writes the DAC offset of the STG for a single channel to permanent EEProm memory

void ErasePermanentDacOffset (LIH30_DAC_Channel_EnumNet DacChannel)

Delets the DAC offset of the STG for a single channel in permanent EEProm memory

void SetAudioOutDacParameter (uint32_t Frequency, uint32_t Amplification)

Sets the parameter of the audio DAC output.

 void GetAudioOutDacParameter ([System::Runtime::InteropServices::Out]uint32_t% Frequency, [System::← Runtime::InteropServices::Out]uint32_t% Amplification)

Gets the parameter of the audio DAC output.

Properties

• CStimulusFunctionNet[^] StimulusFunction [get]

Additional Inherited Members

11.39.1 Detailed Description

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

11.39.2 Constructor & Destructor Documentation

```
11.39.2.1 CLIH3DeviceNet() CLIH3DeviceNet ()
```

Initializes a new instance of the CLIH3DeviceNet class.

```
11.39.2.2 ~CLIH3DeviceNet() virtual ~CLIH3DeviceNet () [virtual]
```

```
11.39.2.3 "!CLIH3DeviceNet() !CLIH3DeviceNet ()
```

11.39.3 Member Function Documentation

```
11.39.3.1 DummyCommand() void DummyCommand (
uint32_t dummyParameter)
```

Dummy command to show how to use the DLL.

Parameters

dummyParameter | parameter to send to the device

```
11.39.3.2 EnableUserTrigger() void EnableUserTrigger ( bool enable )
```

Enables the User Trigger

Parameters

enable Enable

Delets the ADC offset of the DACQ for a single channel in permanent EEProm memory

Parameters

AdcChannel The ADC of	channel
-----------------------	---------

Delets the DAC offset of the STG for a single channel in permanent EEProm memory

Parameters

DacChannel	The DAC channel
DacChannei	The DAC channel

Gets the ADC offset of the DACQ for a single channel

Parameters

AdcChannel	The ADC channel

Returns

The offset for the given channel number

Gets the parameter of the audio DAC output.

Frequency	Frequency(1 - 25000 Hz)
Amplification	Amplification(0 - 0xFFFF)

```
11.39.3.7 GetDacIdleValue() int32_t GetDacIdleValue ( uint32_t DacChannel )
```

Gets the DAC Idle value

Parameters

```
DacChannel The DAC channel
```

Returns

The idle value

Gets the offset of a DAC channel.

Parameters

```
DacChannel The DAC channel
```

Returns

The offset for the given channel number

```
11.39.3.9 GetDacqRunStatus() uint16_t GetDacqRunStatus ( )
```

Gets the data acquisition running status

Returns

The status (1: running / 0: stopped)

```
11.39.3.10 GetDacUseldleValue() bool GetDacUseIdleValue ( uint32_t DacChannel )
```

Gets if the DAC Idle value is used after stimulation

e DAC channel

Returns

Use idle value

11.39.3.11 GetDigInState() uint16_t GetDigInState ()

Reads from the LIH30 digital input

Returns

The bit mask defining the digital input state

Reads the requested amount of EEprom byte from the EPC10 EEPROM

Parameters

EEpromStartAddress	start address of memory area to read from
EEpromData_Length	The maximal length of EEpromData.
epc10bus	The EPC10 bus

Returns

pointer to internal memory for the requested amount of data

11.39.3.13 GetSampleInterval() uint32_t GetSampleInterval ()

Gets the Sample Interval for the DACQ and Stimulation

Returns

Sample Interval configured on the device

11.39.3.14 IsUserTriggerEnabled() bool IsUserTriggerEnabled ()

Is the User Trigger enabled

Returns

Enabled

Gets the clipping information

Parameters

```
epc10bus The EPC10 bus
```

Returns

The clipping value

Send command to the EPC10

Parameters

epc10bus	The EPC10 bus
Command	The command

Sets the ADC offset of the DACQ for a single channel

AdcChannel	The ADC channel
Offset	The offset for the given channel number

Writes the ADC offset of the DACQ for a single channel to permanent EEProm memory

Parameters

AdcChannel	The ADC channel
riaconanion	THO ADO GHAINIO

Sets the parameter of the audio DAC output.

Parameters

Frequency	Frequency(1 - 25000 Hz)
Amplification	Amplification(0 - 0xFFFF)

Sets the DAC Idle value

Parameters

DacChannel	The DAC channel
IdleValue	The idle value

Sets the offset of a DAC channel.

DacChannel	The DAC channel
Offset	The offset for the given channel number

Writes the DAC offset of the STG for a single channel to permanent EEProm memory

Parameters

DacChannel	The DAC channel
Daoonanion	The Brite dilatine

Sets if the DAC Idle value is used after stimulation

Parameters

DacChannel	The DAC channel
Useldle	Use idle value

```
11.39.3.24 SetDigOutState() void SetDigOutState ( uint16_t DigOutState )
```

Writes to the LIH30 digital output

Parameters

DigOutState	The bit mask defining the digital output state
-------------	--

Writes into EEprom on the EPC10 EEPROM

EEpromStartAddress	start address of memory area to write to
EEpromData	pointer to internal memory for the supported amount of data
epc10bus	The EPC10 bus

```
11.39.3.26 SetSampleInterval() void SetSampleInterval ( uint32_t SampleInterval)
```

Sets the Sample Interval for the DACQ and Stimulation

Parameters


	SampleInterval	between the samples, Sample interval is available from 1 to 4194303
--	----------------	---

11.39.4 Property Documentation

11.39.4.1 StimulusFunction CStimulusFunctionNet^ StimulusFunction [get]

11.40 CMcsBus_AxisParametersNet Class Reference

Inheritance diagram for CMcsBus AxisParametersNet:

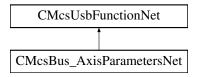

Public Member Functions

- CMcsBus_AxisParametersNet (CMcsUsbNet[^] device)
- ~CMcsBus_AxisParametersNet (void)
- void SetAxisParametersEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, unsigned int parameter)
- void SetAxisParametersEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, int parameter)
- unsigned int GetAxisParametersUnsignedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)
- int GetAxisParametersSignedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)


Additional Inherited Members

11.40.1 Constructor & Destructor Documentation

```
11.40.1.1 CMcsBus_AxisParametersNet() CMcsBus_AxisParametersNet (
             CMcsUsbNet^ device )
11.40.1.2 ~ CMcsBus_AxisParametersNet() ~ CMcsBus_AxisParametersNet (
             void )
11.40.2 Member Function Documentation
11.40.2.1 GetAxisParametersSignedEeprom() int GetAxisParametersSignedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short index)
11.40.2.2 GetAxisParametersUnsignedEeprom() unsigned int GetAxisParametersUnsignedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short index )
11.40.2.3 SetAxisParametersEeprom() [1/2] void SetAxisParametersEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short index,
             int parameter )
11.40.2.4 SetAxisParametersEeprom() [2/2] void SetAxisParametersEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short index,
             unsigned int parameter )
```

11.41 CMcsBus_ExtensionNet Class Reference

Inheritance diagram for CMcsBus_ExtensionNet:

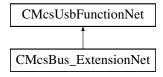

Public Member Functions

- CMcsBus_ExtensionNet (CMcsUsbNet[^] device)
- ∼CMcsBus_ExtensionNet (void)
- · void SetLEDSwitch (unsigned char busnumber, unsigned char busaddress, unsigned short LEDSwitch)
- unsigned short GetLEDSwitch (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.41.1 Constructor & Destructor Documentation


11.41.2 Member Function Documentation

```
11.41.2.1 GetLEDSwitch() unsigned short GetLEDSwitch (
unsigned char busnumber,
unsigned char busaddress)
```

```
11.41.2.2 SetLEDSwitch() void SetLEDSwitch (
unsigned char busnumber,
unsigned char busaddress,
unsigned short LEDSwitch)
```

11.42 CMcsBus_FYIExtensionNet Class Reference

Inheritance diagram for CMcsBus_FYIExtensionNet:

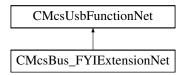

Public Member Functions

- CMcsBus_FYIExtensionNet (CMcsUsbNet[^] device)
- ~CMcsBus_FYIExtensionNet (void)
- void SetValves (unsigned char busnumber, unsigned char busaddress, unsigned int states)
- unsigned int GetValves (unsigned char busnumber, unsigned char busaddress)
- void SetDIO (unsigned char busnumber, unsigned char busaddress, unsigned short io)
- unsigned short GetDIO (unsigned char busnumber, unsigned char busaddress)
- void SetSingleHeater (unsigned char busnumber, unsigned char busaddress, short index, unsigned short power)
- unsigned short GetSingleHeater (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members

11.42.1 Constructor & Destructor Documentation

```
11.42.2.2 GetSingleHeater() unsigned short GetSingleHeater (
    unsigned char busnumber,
    unsigned char busaddress,
    short index )
```

```
11.42.2.3 GetValves() unsigned int GetValves (
unsigned char busnumber,
unsigned char busaddress)
```


```
11.42.2.4 SetDIO() void SetDIO (
          unsigned char busnumber,
          unsigned char busaddress,
          unsigned short io )

11.42.2.5 SetSingleHeater() void SetSingleHeater (
          unsigned char busnumber,
          unsigned char busaddress,
          short index,
          unsigned short power )

11.42.2.6 SetValves() void SetValves (
          unsigned char busnumber,
          unsigned char busaddress,
          unsigned int states )
```

11.43 CMcsBus MotorControlNet Class Reference

Inheritance diagram for CMcsBus MotorControlNet:

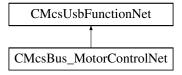

- CMcsBus_MotorControlNet (CMcsUsbNet[^] device)
- ~CMcsBus_MotorControlNet (void)
- void SetMCScalingFactorEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int GetMCScalingFactorEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCScalingFactor (unsigned char busnumber, unsigned char busnumber, under bu
- int GetMCScalingFactor (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short GetMCMaxSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short GetMCMaxSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravelEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int GetMCMaxTravelEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravel (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)

- int GetMCMaxTravel (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- · void SetMCCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCMaxCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCMaxCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRegulatorGainEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short GetMCRegulatorGainEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRegulatorGain (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short GetMCRegulatorGain (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCMaxAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCMaxAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short GetMCStandbyCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short GetMCStandbyCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyTimeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short GetMCStandbyTimeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyTime (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short GetMCStandbyTime (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCBreakCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCBreakCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCBreakCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCBreakCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCConfigEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short GetMCConfigEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCConfig (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short confia)
- unsigned short GetMCConfig (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)

- unsigned short GetMCSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short GetMCSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReferenceCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCReferenceCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReferenceCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCReferenceCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentModeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet GetMCCurrentModeEeprom (unsigned char busnumber, unsigned char axis)
- void SetMCCurrentMode (unsigned char busnumber, unsigned char busaddress, unsigned char axis, Robo
 —
 CurrentModeEnumNet mode)
- RoboCurrentModeEnumNetGetMCCurrentMode (unsigned char busnumber, unsigned char busn
- void SetMCAxisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short revision)
- unsigned short GetMCAxisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedUnitEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int32_t speedunit)
- int32 t GetMCSpeedUnitEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char axis, bool OnOff_status)
- bool GetMCOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short GetMCSpeedShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrentShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravelShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int GetMCMaxTravelShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int GetMCCurrentPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void SetMCNewPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int GetMCNewPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- short GetMCCurrentSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void StartMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRotation (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char onoff)
- unsigned short GetMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReference (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char switch_enable, unsigned char switch_polarity)
- unsigned char GetMCReference (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned char% switch_port)
- void StopMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentModeShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet GetMCCurrentModeShortCommand (unsigned char busnumber, unsigned char busnumber, unsigned char axis)
- unsigned short GetMCPhase (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short GetMCPhaseOffset (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetSubChannel (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short subchannel)
- unsigned short GetSubChannel (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

Additional Inherited Members

11.43.1 Constructor & Destructor Documentation

11.43.2 Member Function Documentation

```
11.43.2.1 GetMCAcceleration() unsigned short GetMCAcceleration (
    unsigned char busnumber,
    unsigned char busaddress,
    unsigned char axis)
```

```
\textbf{11.43.2.2} \quad \textbf{GetMCAccelerationEeprom()} \quad \textbf{unsigned short GetMCAccelerationEeprom ()}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
\textbf{11.43.2.3} \quad \textbf{GetMCAccelerationShortCommand()} \quad \textbf{unsigned short GetMCAccelerationShortCommand ()} \\
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
\textbf{11.43.2.4} \quad \textbf{GetMCAxisRevisionEeprom()} \quad \textbf{unsigned short GetMCAxisRevisionEeprom ()}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
11.43.2.5 GetMCBreakCurrent() short GetMCBreakCurrent (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
11.43.2.6 GetMCBreakCurrentEeprom() short GetMCBreakCurrentEeprom (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
11.43.2.7 GetMCConfig() unsigned short GetMCConfig (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
11.43.2.8 GetMCConfigEeprom() unsigned short GetMCConfigEeprom (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
```

```
11.43.2.9 GetMCCurrent() short GetMCCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.10 GetMCCurrentEeprom() short GetMCCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.11 GetMCCurrentMode() RoboCurrentModeEnumNet GetMCCurrentMode (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.12 GetMCCurrentModeEeprom() RoboCurrentModeEnumNet GetMCCurrentModeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.13 GetMCCurrentModeShortCommand() RoboCurrentModeEnumNet GetMCCurrentModeShort←
Command (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.14 GetMCCurrentPosition() int GetMCCurrentPosition (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.15 GetMCCurrentShortCommand() short GetMCCurrentShortCommand (
            unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.43.2.16 GetMCCurrentSpeed() short GetMCCurrentSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.17 GetMCMaxAcceleration() unsigned short GetMCMaxAcceleration (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.18 GetMCMaxAccelerationEeprom() unsigned short GetMCMaxAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.19 GetMCMaxCurrent() short GetMCMaxCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.20 GetMCMaxCurrentEeprom() short GetMCMaxCurrentEeprom (
            unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.21 GetMCMaxSpeed() unsigned short GetMCMaxSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.22 GetMCMaxSpeedEeprom() unsigned short GetMCMaxSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.43.2.23 GetMCMaxTravel() int GetMCMaxTravel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.24 GetMCMaxTravelEeprom() int GetMCMaxTravelEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.25 GetMCMaxTravelShortCommand() int GetMCMaxTravelShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.26 GetMCMovement() unsigned short GetMCMovement (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.27 GetMCNewPosition() int GetMCNewPosition (
            unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.28 GetMCOutputOnOff() bool GetMCOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.29 GetMCPhase() unsigned short GetMCPhase (
            unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.43.2.30 GetMCPhaseOffset() unsigned short GetMCPhaseOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.31 GetMCReference() unsigned char GetMCReference (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned char% switch_port)
11.43.2.32 GetMCReferenceCurrent() short GetMCReferenceCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.33 GetMCReferenceCurrentEeprom() short GetMCReferenceCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.34 GetMCRegulatorGain() short GetMCRegulatorGain (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.35 GetMCRegulatorGainEeprom() short GetMCRegulatorGainEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.36 GetMCScalingFactor() int GetMCScalingFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.43.2.37 GetMCScalingFactorEeprom() int GetMCScalingFactorEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.38 GetMCSpeed() short GetMCSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.39 GetMCSpeedEeprom() unsigned short GetMCSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.40 GetMCSpeedShortCommand() short GetMCSpeedShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.41 GetMCSpeedUnitEeprom() int32_t GetMCSpeedUnitEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.42 GetMCStandbyCurrent() short GetMCStandbyCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.43 GetMCStandbyCurrentEeprom() short GetMCStandbyCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.43.2.44 GetMCStandbyTime() short GetMCStandbyTime (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.45 GetMCStandbyTimeEeprom() short GetMCStandbyTimeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.46 GetSubChannel() unsigned short GetSubChannel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.43.2.47 SetMCAcceleration() void SetMCAcceleration (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.43.2.48 SetMCAccelerationEeprom() void SetMCAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.43.2.49 SetMCAccelerationShortCommand() void SetMCAccelerationShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.43.2.50 SetMCAxisRevisionEeprom() void SetMCAxisRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short revision )
```

```
11.43.2.51 SetMCBreakCurrent() void SetMCBreakCurrent (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short current )
\textbf{11.43.2.52} \quad \textbf{SetMCBreakCurrentEeprom()} \quad \texttt{void SetMCBreakCurrentEeprom} \quad \textbf{(}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short current )
11.43.2.53 SetMCConfig() void SetMCConfig (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              unsigned short config )
11.43.2.54 SetMCConfigEeprom() void SetMCConfigEeprom (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              unsigned short config )
11.43.2.55 SetMCCurrent() void SetMCCurrent (
             unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short current )
11.43.2.56 SetMCCurrentEeprom() void SetMCCurrentEeprom (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short current )
```

```
11.43.2.57 SetMCCurrentMode() void SetMCCurrentMode (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              {\tt RoboCurrentModeEnumNet}\ \textit{mode}\ )
\textbf{11.43.2.58} \quad \textbf{SetMCCurrentModeEeprom()} \quad \texttt{void SetMCCurrentModeEeprom} \quad \textbf{(}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              RoboCurrentModeEnumNet mode )
11.43.2.59 SetMCCurrentModeShortCommand() void SetMCCurrentModeShortCommand (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              RoboCurrentModeEnumNet mode )
11.43.2.60 SetMCCurrentPosition() void SetMCCurrentPosition (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              int position )
11.43.2.61 SetMCCurrentShortCommand() void SetMCCurrentShortCommand (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short current )
11.43.2.62 SetMCMaxAcceleration() void SetMCMaxAcceleration (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              unsigned short acceleration )
```

```
\textbf{11.43.2.63} \quad \textbf{SetMCMaxAccelerationEeprom()} \quad \texttt{void SetMCMaxAccelerationEeprom} \quad \textbf{(}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.43.2.64 SetMCMaxCurrent() void SetMCMaxCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.43.2.65 SetMCMaxCurrentEeprom() void SetMCMaxCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.43.2.66 SetMCMaxSpeed() void SetMCMaxSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
11.43.2.67 SetMCMaxSpeedEeprom() void SetMCMaxSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
11.43.2.68 SetMCMaxTravel() void SetMCMaxTravel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
```

```
11.43.2.69 SetMCMaxTravelEeprom() void SetMCMaxTravelEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
11.43.2.70 SetMCMaxTravelShortCommand() void SetMCMaxTravelShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
11.43.2.71 SetMCNewPosition() void SetMCNewPosition (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int position )
11.43.2.72 SetMCOutputOnOff() void SetMCOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             bool OnOff_status )
11.43.2.73 SetMCReference() void SetMCReference (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned char switch_enable,
             unsigned char switch_polarity )
11.43.2.74 SetMCReferenceCurrent() void SetMCReferenceCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
```

```
11.43.2.75 SetMCReferenceCurrentEeprom() void SetMCReferenceCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
\textbf{11.43.2.76} \quad \textbf{SetMCRegulatorGain()} \quad \texttt{void SetMCRegulatorGain ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short gain )
11.43.2.77 SetMCRegulatorGainEeprom() void SetMCRegulatorGainEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short gain )
11.43.2.78 SetMCRotation() void SetMCRotation (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned char onoff )
11.43.2.79 SetMCScalingFactor() void SetMCScalingFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int factor )
11.43.2.80 SetMCScalingFactorEeprom() void SetMCScalingFactorEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int factor )
```


```
11.43.2.81 SetMCSpeed() void SetMCSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short speed )
11.43.2.82 SetMCSpeedEeprom() void SetMCSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
11.43.2.83 SetMCSpeedShortCommand() void SetMCSpeedShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short speed )
11.43.2.84 SetMCSpeedUnitEeprom() void SetMCSpeedUnitEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int32_t speedunit )
11.43.2.85 SetMCStandbyCurrent() void SetMCStandbyCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short percent )
11.43.2.86 SetMCStandbyCurrentEeprom() void SetMCStandbyCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short percent )
```

```
11.43.2.87 SetMCStandbyTime() void SetMCStandbyTime (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short t )
\textbf{11.43.2.88} \quad \textbf{SetMCS} \\ \textbf{tandbyTimeEeprom()} \quad \texttt{void SetMCS} \\ \textbf{tandbyTimeEeprom ()}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              short t )
11.43.2.89 SetSubChannel() void SetSubChannel (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              unsigned short subchannel )
11.43.2.90 StartMCMovement() void StartMCMovement (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis )
11.43.2.91 StopMCMovement() void StopMCMovement (
              unsigned char busnumber,
              unsigned char busaddress,
```

11.44 CMcsBus_SensorNet Class Reference

unsigned char axis)

Inheritance diagram for CMcsBus_SensorNet:

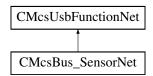

- CMcsBus SensorNet (CMcsUsbNet[∧] device)
- ~CMcsBus_SensorNet (void)
- · void SetMinimalThreshold (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short GetMinimalThreshold (unsigned char busnumber, unsigned char busaddress)
- void SetDetectionThreshold (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short GetDetectionThreshold (unsigned char busnumber, unsigned char busaddress)
- void SetLatency (unsigned char busnumber, unsigned char busaddress, unsigned short latency)
- unsigned short GetLatency (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetBubbleStatus (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetLatencyCounter (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetDetectorValue (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ GetPressure (unsigned char busnumber, unsigned char busaddress, int n)
- int GetPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRegulatorOnOff (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned char onoff)
- unsigned char GetRegulatorOnOff (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetSollPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index, int pressure)
- int GetSollPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRegulatorFactor (unsigned char busnumber, unsigned char busaddress, unsigned short index, int factor)
- int GetRegulatorFactor (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetPressureOffset (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- array< unsigned short > ^ GetPressureOffset (unsigned char busnumber, unsigned char busaddress)
- int GetPressureOffset (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- unsigned int GetRegulatorStatus (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRotatePump (unsigned char busnumber, unsigned char busaddress, unsigned short index, short speed)
- short GetRotatePump (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetMovePump (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned short speed, int position)

- array< int > ^ Get4ADC (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ Get4ADCAverage (unsigned char busnumber, unsigned char busaddress)
- void Set4DAC (unsigned char busnumber, unsigned char busaddress, array< unsigned short >^ dac)
- array< unsigned short > ^ Get4DAC (unsigned char busnumber, unsigned char busaddress)
- void Set4ADCMode (unsigned char busnumber, unsigned char busaddress, PatchServAdcModeEnumNet mode)
- PatchServAdcModeEnumNet Get4ADCMode (unsigned char busnumber, unsigned char busaddress)
- void Set4ADCCatchampAverageShift (unsigned char busnumber, unsigned char busaddress, unsigned int shift)
- unsigned int Get4ADCCatchampAverageShift (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ Get2AnalogInput (unsigned char busnumber, unsigned char busaddress)
- unsigned short Get2DigitalInput (unsigned char busnumber, unsigned char busaddress)

- array< unsigned short > ^ GetADCs (unsigned char busnumber, unsigned char busaddress, int n)
- array< unsigned short > ^ GetADCsLoop (unsigned char busnumber, unsigned char busaddress, int n)
- void SetPiezoState (unsigned char busnumber, unsigned char busaddress, int state)
- void SetDACs (unsigned char busnumber, unsigned char busaddress, unsigned short index, array< unsigned short >^ dac_times_voltages)
- array< unsigned short > ^ GetDACs (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- · void SetSamplePeriode (unsigned char busnumber, unsigned char busaddress, unsigned short periode)
- unsigned short GetSamplePeriode (unsigned char busnumber, unsigned char busaddress)
- void StartSync (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetSyncState (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacAmplitude (unsigned char busnumber, unsigned char busaddress, unsigned short dacAmplitude)
- unsigned short CatchAmpGetDacAmplitude (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacOffset (unsigned char busnumber, unsigned char busaddress, short dacOffset)
- short CatchAmpGetDacOffset (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcMean (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcValue (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcValueH (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcValueL (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetPwmEnable (unsigned char busnumber, unsigned char busaddress, bool pwmEnable)
- bool CatchAmpGetPwmEnable (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacEnable (unsigned char busnumber, unsigned char busaddress, bool dacEnable)
- bool CatchAmpGetDacEnable (unsigned char busnumber, unsigned char busaddress)
- int TactSwitchGetState (unsigned char busnumber, unsigned char busaddress)
- void TactSwitchSetDisplay (unsigned char busnumber, unsigned char busaddress, int Melody)

Additional Inherited Members

11.44.1 Constructor & Destructor Documentation

11.44.2 Member Function Documentation

```
11.44.2.1 CatchAmpGetAdcMean() int CatchAmpGetAdcMean (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.2 CatchAmpGetAdcValue() int CatchAmpGetAdcValue (
             unsigned char busnumber,
             unsigned char busaddress )
\textbf{11.44.2.3} \quad \textbf{CatchAmpGetAdcValueH()} \quad \texttt{int CatchAmpGetAdcValueH ()}
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.4 CatchAmpGetAdcValueL() int CatchAmpGetAdcValueL (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.5 CatchAmpGetDacAmplitude() unsigned short CatchAmpGetDacAmplitude (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.6 CatchAmpGetDacEnable() bool CatchAmpGetDacEnable (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.7 CatchAmpGetDacOffset() short CatchAmpGetDacOffset (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.8 CatchAmpGetPwmEnable() bool CatchAmpGetPwmEnable (
             unsigned char busnumber,
             unsigned char busaddress )
```

```
11.44.2.9 CatchAmpSetDacAmplitude() void CatchAmpSetDacAmplitude (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short dacAmplitude )
11.44.2.10 CatchAmpSetDacEnable() void CatchAmpSetDacEnable (
             unsigned char busnumber,
             unsigned char busaddress,
             bool dacEnable )
11.44.2.11 CatchAmpSetDacOffset() void CatchAmpSetDacOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             short dacOffset )
11.44.2.12 CatchAmpSetPwmEnable() void CatchAmpSetPwmEnable (
             unsigned char busnumber,
             unsigned char busaddress,
             bool pwmEnable )
11.44.2.13 Get2AnalogInput() array<unsigned short> ^ Get2AnalogInput (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.14 Get2DigitalInput() unsigned short Get2DigitalInput (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.15 Get4ADC() array<int> ^{\land} Get4ADC (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.16 Get4ADCAverage() array<int> ^ Get4ADCAverage (
             unsigned char busnumber,
             unsigned char busaddress )
```

```
\textbf{11.44.2.17} \quad \textbf{Get4ADCCatchampAverageShift()} \quad \textbf{unsigned int Get4ADCCatchampAverageShift ()}
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.18 Get4ADCMode() PatchServAdcModeEnumNet Get4ADCMode (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.19 Get4DAC() array<unsigned short> ^{\land} Get4DAC (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.20 GetADCs() array<unsigned short> ^ GetADCs (
             unsigned char busnumber,
             unsigned char busaddress,
             int n)
11.44.2.21 GetADCsLoop() array<unsigned short> ^{\land} GetADCsLoop (
             unsigned char busnumber,
             unsigned char busaddress,
             int n)
11.44.2.22 GetBubbleStatus() unsigned short GetBubbleStatus (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.23 GetDACs() array<unsigned short> ^{\land} GetDACs (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
11.44.2.24 GetDetectionThreshold() unsigned short GetDetectionThreshold (
             unsigned char busnumber,
             unsigned char busaddress )
```

```
{\bf 11.44.2.25} \quad {\bf GetDetectorValue()} \quad {\tt unsigned short GetDetectorValue ()}
              unsigned char busnumber,
              unsigned char busaddress )
11.44.2.26 GetLatency() unsigned short GetLatency (
              unsigned char busnumber,
              unsigned char busaddress )
11.44.2.27 GetLatencyCounter() unsigned short GetLatencyCounter (
              unsigned char busnumber,
              unsigned char busaddress )
\textbf{11.44.2.28} \quad \textbf{GetMinimalThreshold()} \quad \texttt{unsigned short GetMinimalThreshold ()}
              unsigned char busnumber,
              unsigned char busaddress )
11.44.2.29 GetMovePump() void GetMovePump (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short index,
              [System::Runtime::InteropServices::Out] unsigned short% speed,
              [System::Runtime::InteropServices::Out] int% position )
11.44.2.30 GetPiezoState() void GetPiezoState (
              unsigned char busnumber,
              unsigned char busaddress,
              [System::Runtime::InteropServices::Out] int% state,
              [System::Runtime::InteropServices::Out] int% reason )
11.44.2.31 GetPressure() [1/2] array<int> ^ GetPressure (
              unsigned char busnumber,
              unsigned char busaddress,
              int n)
```

```
11.44.2.32 GetPressure() [2/2] int GetPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.44.2.33 GetPressureOffset() [1/2] array<unsigned short> ^ GetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress )
11.44.2.34 GetPressureOffset() [2/2] int GetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.44.2.35 GetRegulationTimeouts() void GetRegulationTimeouts (
             unsigned char busnumber,
             unsigned char busaddress,
             [System::Runtime::InteropServices::Out] unsigned short% MaxSpeedWait,
             [System::Runtime::InteropServices::Out] unsigned short% MaxSignChange )
11.44.2.36 GetRegulatorFactor() int GetRegulatorFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.44.2.37 GetRegulatorOnOff() unsigned char GetRegulatorOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.44.2.38 GetRegulatorStatus() unsigned int GetRegulatorStatus (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
```

```
\textbf{11.44.2.39} \quad \textbf{GetRotatePump()} \quad \texttt{short GetRotatePump} \quad \textbf{(}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short index)
{\bf 11.44.2.40} \quad {\bf GetSamplePeriode()} \quad {\tt unsigned short GetSamplePeriode } \ (
              unsigned char busnumber,
              unsigned char busaddress )
11.44.2.41 GetSollPressure() int GetSollPressure (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short index )
11.44.2.42 GetSyncState() unsigned short GetSyncState (
              unsigned char busnumber,
              unsigned char busaddress )
11.44.2.43 Set4ADCCatchampAverageShift() void Set4ADCCatchampAverageShift (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned int shift )
11.44.2.44 Set4ADCMode() void Set4ADCMode (
              unsigned char busnumber,
              unsigned char busaddress,
              PatchServAdcModeEnumNet mode )
11.44.2.45 Set4DAC() void Set4DAC (
              unsigned char busnumber,
              unsigned char busaddress,
              array< unsigned short >^{\wedge} dac )
```


```
11.44.2.46 SetDACs() void SetDACs (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             array< unsigned short >^{\wedge} dac\_times\_voltages )
11.44.2.47 SetDetectionThreshold() void SetDetectionThreshold (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short threshold )
11.44.2.48 SetLatency() void SetLatency (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short latency )
11.44.2.49 SetMinimalThreshold() void SetMinimalThreshold (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short threshold )
11.44.2.50 SetMovePump() void SetMovePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             unsigned short speed,
             int position )
11.44.2.51 SetPiezoState() void SetPiezoState (
             unsigned char busnumber,
             unsigned char busaddress,
             int state )
11.44.2.52 SetPressureOffset() void SetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
```

```
11.44.2.53 SetRegulationTimeouts() void SetRegulationTimeouts (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short MaxSpeedWait,
             unsigned short MaxSignChange )
11.44.2.54 SetRegulatorFactor() void SetRegulatorFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             int factor )
11.44.2.55 SetRegulatorOnOff() void SetRegulatorOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             unsigned char onoff )
11.44.2.56 SetRotatePump() void SetRotatePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             short speed )
11.44.2.57 SetSamplePeriode() void SetSamplePeriode (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short periode )
11.44.2.58 SetSollPressure() void SetSollPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             int pressure )
11.44.2.59 StartSync() void StartSync (
             unsigned char busnumber,
             unsigned char busaddress )
```

11.45 CMcsBus_TempSensorNet Class Reference

Inheritance diagram for CMcsBus_TempSensorNet:

int Melody)

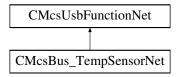

Public Member Functions

- CMcsBus_TempSensorNet (CMcsUsbNet[∧] device)
- ~CMcsBus TempSensorNet (void)
- short GetTemperatur (unsigned char busnumber, unsigned char busaddress)
- short GetTemperatur (unsigned char busnumber, unsigned char busaddress, short index)
- void SetNanoVoltsPerKelvin (unsigned char busnumber, unsigned char busaddress, int nanovoltsperkelvin)
- int GetNanoVoltsPerKelvin (unsigned char busnumber, unsigned char busaddress)
- short GetThermoVoltage (unsigned char busnumber, unsigned char busaddress, short index)
- short GetThermoTemp (unsigned char busnumber, unsigned char busaddress, short index)
- void SetThermoOffset (unsigned char busnumber, unsigned char busaddress, short index, short offset)
- short GetThermoOffset (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members

11.45.1 Constructor & Destructor Documentation


11.45.2 Member Function Documentation

```
11.45.2.1 GetNanoVoltsPerKelvin() int GetNanoVoltsPerKelvin (
             unsigned char busnumber,
             unsigned char busaddress )
11.45.2.2 GetTemperatur() [1/2] short GetTemperatur (
             unsigned char busnumber,
             unsigned char busaddress )
11.45.2.3 GetTemperatur() [2/2] short GetTemperatur (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.45.2.4 GetThermoOffset() short GetThermoOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.45.2.5 GetThermoTemp() short GetThermoTemp (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.45.2.6 GetThermoVoltage() short GetThermoVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.45.2.7 SetNanoVoltsPerKelvin() void SetNanoVoltsPerKelvin (
             unsigned char busnumber,
             unsigned char busaddress,
             int nanovoltsperkelvin )
```

```
11.45.2.8 SetThermoOffset() void SetThermoOffset (
          unsigned char busnumber,
          unsigned char busaddress,
          short index,
          short offset )
```

11.46 CMcsBus_VoltageModeNet Class Reference

Inheritance diagram for CMcsBus_VoltageModeNet:

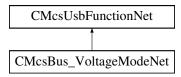

- CMcsBus_VoltageModeNet (CMcsUsbNet[^] device)
- ~CMcsBus VoltageModeNet (void)
- void SetVMMaxPositiveCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxPositiveCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxPositiveCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxNegativeCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxNegativeCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxPositiveVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxPositiveVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxNegativeVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxNegativeVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)

- · void SetVMOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char channel, unsigned short status)
- unsigned short GetVMOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMVoltage (unsigned char busnumber, unsigned char busnumber, under busnumber, unde voltage)
- short GetVMVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)

Additional Inherited Members

11.46.1 Constructor & Destructor Documentation

```
11.46.1.1 CMcsBus_VoltageModeNet() CMcsBus_VoltageModeNet (
            CMcsUsbNet^ device )
11.46.1.2 ~ CMcsBus_VoltageModeNet() ~ CMcsBus_VoltageModeNet (
            void )
11.46.2 Member Function Documentation
```

```
11.46.2.1 GetVMMaxNegativeCurrent() short GetVMMaxNegativeCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
```

```
11.46.2.2 GetVMMaxNegativeCurrentEeprom() short GetVMMaxNegativeCurrentEeprom (
            unsigned char busnumber,
```

```
unsigned char busaddress,
unsigned char channel )
```

unsigned char channel)

11.46.2.3 GetVMMaxNegativeVoltage() short GetVMMaxNegativeVoltage (

```
unsigned char busnumber,
unsigned char busaddress,
unsigned char channel )
```

```
\textbf{11.46.2.4} \quad \textbf{GetVMMaxNegativeVoltageEeprom()} \quad \texttt{short GetVMMaxNegativeVoltageEeprom ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.5 GetVMMaxPositiveCurrent() short GetVMMaxPositiveCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.6 GetVMMaxPositiveCurrentEeprom() short GetVMMaxPositiveCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.7 GetVMMaxPositiveVoltage() short GetVMMaxPositiveVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.8 GetVMMaxPositiveVoltageEeprom() short GetVMMaxPositiveVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.9 GetVMOutputOnOff() unsigned short GetVMOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.46.2.10 GetVMVoltage() short GetVMVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
```


```
11.46.2.11 SetVMMaxNegativeCurrent() void SetVMMaxNegativeCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.46.2.12 SetVMMaxNegativeCurrentEeprom() void SetVMMaxNegativeCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.46.2.13 SetVMMaxNegativeVoltage() void SetVMMaxNegativeVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
11.46.2.14 SetVMMaxNegativeVoltageEeprom() void SetVMMaxNegativeVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
11.46.2.15 SetVMMaxPositiveCurrent() void SetVMMaxPositiveCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.46.2.16 SetVMMaxPositiveCurrentEeprom() void SetVMMaxPositiveCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
```

11.46.2.20 SetVMVoltage() void SetVMVoltage (

```
unsigned char busnumber,
unsigned char busaddress,
unsigned char channel,
short voltage )
```

11.47 CMcsBusNet Class Reference

Inheritance diagram for CMcsBusNet:

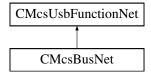

- CMcsBusNet (CMcsUsbNet[^] device)
- virtual ~CMcsBusNet (void)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short value)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, short value)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned int value)

- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, int value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned short% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]short% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned int% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]int% value)
- void SetBusAddressEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short GetBusAddressEeprom (unsigned char busnumber, unsigned char busaddress)
- void SetBusAddress (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short GetBusAddress (unsigned char busnumber, unsigned char busaddress)
- void CMcsBusNet::SetModeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short CMcsBusNet::GetModeEeprom (unsigned char busnumber, unsigned char busaddress)
- void CMcsBusNet::SetMode (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short CMcsBusNet::GetMode (unsigned char busnumber, unsigned char busaddress)
- void SetHWRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short revision)
- unsigned short GetHWRevisionEeprom (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.47.1 Constructor & Destructor Documentation

11.47.2 Member Function Documentation

```
\textbf{11.47.2.2} \quad \textbf{CMcsBusNet::} \textbf{GetModeEeprom()} \quad \textbf{unsigned short CMcsBusNet::} \textbf{GetModeEeprom ()}
              unsigned char busnumber,
              unsigned char busaddress )
11.47.2.3 CMcsBusNet::SetMode() void CMcsBusNet::SetMode (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short mode )
11.47.2.4 CMcsBusNet::SetModeEeprom() void CMcsBusNet::SetModeEeprom (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short mode )
11.47.2.5 GetBusAddress() unsigned short GetBusAddress (
              unsigned char busnumber,
              unsigned char busaddress )
\textbf{11.47.2.6} \quad \textbf{GetBusAddressEeprom()} \quad \texttt{unsigned short GetBusAddressEeprom ()}
              unsigned char busnumber,
              unsigned char busaddress )
11.47.2.7 GetCommand() [1/4] void GetCommand (
              unsigned char command,
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              [System::Runtime::InteropServices::Out] int% value )
11.47.2.8 GetCommand() [2/4] void GetCommand (
              unsigned char command,
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned char axis,
              [System::Runtime::InteropServices::Out] short% value )
```

```
11.47.2.9 GetCommand() [3/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned int% value )
11.47.2.10 GetCommand() [4/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned short% value )
11.47.2.11 GetHWRevisionEeprom() unsigned short GetHWRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress )
11.47.2.12 SetBusAddress() void SetBusAddress (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short newaddress )
11.47.2.13 SetBusAddressEeprom() void SetBusAddressEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short newaddress )
11.47.2.14 SetCommand() [1/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int value )
11.47.2.15 SetCommand() [2/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short value )
```

11.47.2.16 SetCommand() [3/4] void SetCommand (

```
unsigned char command,
unsigned char busnumber,
unsigned char busaddress,
unsigned char axis,
unsigned int value)
```

11.47.2.17 SetCommand() [4/4] void SetCommand (

```
unsigned char command,
unsigned char busnumber,
unsigned char busaddress,
unsigned char axis,
unsigned short value)
```


11.47.2.18 SetHWRevisionEeprom() void SetHWRevisionEeprom (

```
unsigned char busaddress,
unsigned char busaddress,
unsigned short revision )
```

11.48 CMcsUsbDacqNet Class Reference

Base class for data acquisition devices.

Inheritance diagram for CMcsUsbDacqNet:

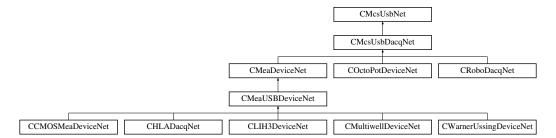

Classes

• class CHWInfo

Class to provide hardware information about the device.

Public Member Functions

- CMcsUsbDacqNet ()
- ∼CMcsUsbDacqNet ()
- virtual uint32_t GetVoltageRangeIndex (unsigned int virtualDevice)
- virtual void SetVoltageRangeByIndex (int32 t voltageRangeIndex, unsigned int virtualDevice)

Sets the voltage range on devices which support multiple voltage ranges.

virtual void SetVoltageRangeInMicroVolt (int32 t voltageRange, unsigned int virtualDevice)

Sets the voltage range on devices which support multiple voltage ranges.

virtual int32 t GetVoltageRangeInMicroVolt (unsigned int virtualDevice)

Gets the currently selected voltage range on devices which support multiple voltage ranges.

virtual int32_t GetVoltageRangeInMilliVolt ()

Gets the currently selected voltage range on devices which support multiple voltage ranges.

virtual void SetDataMode (DataModeEnumNet dataMode, unsigned int virtualDevice)

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

virtual DataModeEnumNet GetDataMode (unsigned int virtualDevice)

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSourceEnumNet source, int bitnumber offset)

Sets the function/source of an digital output bit.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, W2100DigitalSourceEnumNet source, int bitnumber offset)

Sets the function/source of an digital output bit.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, SCUDigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MEA2100_256DigitalSource
 — EnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

• template<typename digitalsourceenum >

void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSource< digital-sourceenum >^ source, int bitnumber_offset)

Sets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::Interop
 — Services::Out]DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber
 — offset)

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::
 — InteropServices::Out]W2100DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::
 — InteropServices::Out]SCUDigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::Interop
 — Services::Out]MEA2100_256DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int%
 bitnumber offset)

Gets the function/source of an digital output bit.

template<typename digitalsourceenum >

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::Interop ← Services::Out]DigitalSource < digitalsourceenum > ^% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

virtual AdapterTypeEnumNet GetAdapterType ()

Gets the adapter which is connected to the MEA2100 device.

virtual MeaLayoutEnumNet GetMeaLayout ()

Gets the MEA layout which is connected to the MEA2100 device.

virtual uint32_t GetAdcDataFormat (uint32_t virtualDevice)

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits,

- virtual uint32_t GetResolutionPerDigit (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% res, [System::Runtime::InteropServices::Out] int% resUnit)
- virtual uint32_t GetHardwareMinRange (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t GetHardwareMaxRange (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t GetDataFormat (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::

 Runtime::InteropServices::Out] int% numberOfBits)
- virtual uint32_t GetNumberOfDataBits (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)

Get the real number of data bits.

virtual void SetSamplerate (int32 t rate, unsigned int oversample, unsigned int virtualDevice)

Sets the sampling frequency of the device.

virtual int32 t GetSamplerate (unsigned int virtualDevice)

Gets the sampling frequency of the device.

virtual uint32_t GetMaxSamplingFrequency (int virtualDevice)

Gets the maximal sampling frequency of the device.

virtual uint32 t GetMinSamplingFrequencyStepsize ()

Gets the minimal sampling frequency step size increment value of the device.

virtual int32_t GetChannelsInBlock (unsigned int virtualDevice)

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

- virtual uint32_t GetChannelLayout ([System::Runtime::InteropServices::Out]int% AnalogChannels, [System::Runtime::InteropServices::Out]int% DigitalChannels, [System::Runtime::InteropServices::Out]int% ChecksumChannels, [System::Runtime::InteropServices::Out]int% TimestampChannels, [System::← Runtime::InteropServices::Out]int% ChannelsInBlock, unsigned int virtualDevice)
- virtual void SendStartDacq ()

Start sampling.

virtual void SendStartDacq (int VirtualDacqMap)

Start sampling.

virtual void SendStartStgAndDacq (uint32_t trigger_map, int VirtualDacqMap)

Start sampling together with the STG.

virtual void SendStopDacq ()

Stop sampling.

virtual void SendStopDacq (int VirtualDacqMap)

Stop sampling.

Parameters

VirtualDacqMap

virtual void SendStopStgAndDacq (uint32_t trigger_map, int VirtualDacqMap)

Stop sampling together with the STG.

• virtual void SendStopStgAndDacqWithOptions (uint32_t trigger_map, int VirtualDacqMap, int options)

Stop sampling together with the STG and options.

virtual void StartLoop ()

Start the data acquisition thread.

virtual void StartLoop (int32 t timeout)

Start the data acquisition thread.

virtual void StartLoop (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_← t packetsInUrb)

Start the data acquisition thread.

virtual void StartLoop (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread.

- virtual void StopLoop ()
- · virtual void ClearBuffers ()
- virtual void StartDacq ()

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout)

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_←
t packetsInUrb)

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread and sampling.

virtual void StopDacq ()

Stop the data acquisition thread and sampling.

virtual void StopDacq (uint32_t virtualDevice)

Stop the data acquisition thread and sampling.

- virtual uint32_t SetPoti (uint32_t channel, uint32_t value, bool write_nvram)
- virtual uint32_t GetPoti (uint32_t channel, [System::Runtime::InteropServices::Out]uint32_t% value)
- virtual CFilterPropertyNet ^ GetFilterProperty (DacqGroupChannelEnumNet GroupID, unsigned int index)
- virtual array< CFilterPropertyNet[^]> [^] CMcsUsbDacqNet::GetFilterProperties (DacqGroupChannelEnum← Net GroupID)
- int GetChannelDataFillSize ()
- virtual void SetSelectedChannels (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void SetSelectedChannels (int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedChannels (array < bool > ^ selectedChannels, int queuesize, int threshold, Sample ← SizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void SetSelectedChannels (array < bool > ^ selectedChannels, int queuesize, int threshold, Sample ← SizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedData (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void SetSelectedData (int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedData (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void SetSelectedData (array < bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read← FramesDict... with handle = 0 to read the data.

- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSizeNet sampleDstSize)
- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, array< bool >[^] selected←
 Channels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read←FramesDict... with handle = 0 to read the data.

- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, array< bool >[∧] selected←
 Channels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSizeNet
- virtual void SetSelectedChannelsQueue (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_\cup ReadFramesDict... with handle = 0 to read the data.

- virtual void SetSelectedChannelsQueue (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedChannelsQueue (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

- virtual void SetSelectedChannelsQueue (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual uint32_t ChannelBlock_AvailFrames (int handle)

Get the number of sample frames already available in the FIFO.

- virtual uint32_t ChannelBlock_AvailFrames (int handle, int queue)

Read data from a FIFO queue in uint16_t data format

virtual uint32_t ChannelBlock_ReadFramesUI16 (int handle, array< uint16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format

virtual array< int16_t > ^ ChannelBlock_ReadFramesI16 (int handle, int frames, [System::Runtime::
 — InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format

virtual uint32_t ChannelBlock_ReadFramesI16 (int handle, array< int16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format

Read data from a FIFO queue in uint32_t data format

virtual uint32_t ChannelBlock_ReadFramesUl32 (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32 t data format

virtual array< int32_t > ^ ChannelBlock_ReadFramesI32 (int handle, int frames, [System::Runtime::
 — InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32 t data format

virtual uint32_t ChannelBlock_ReadFramesl32 (int handle, array< int32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint32_t data format

• virtual array< array< uint16_t >^> ^ ChannelBlock_ReadAsFrameArrayUI16 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< uint16_t >^> ^ ChannelBlock_ReadAsFrameArrayUI16 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

virtual array< array< int16_t >^> ^ ChannelBlock_ReadAsFrameArrayl16 (int handle, int frames, [System← ::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< int16_t >^> ^ ChannelBlock_ReadAsFrameArrayI16 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< uint32_t >^> ^ ChannelBlock_ReadAsFrameArrayUI32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< uint32_t >^> ^ ChannelBlock_ReadAsFrameArrayUl32 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16 t data frame arrays

virtual array< array< int32_t >^> ^ ChannelBlock_ReadAsFrameArrayl32 (int handle, int frames, [System ← ::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< int32_t >^> ^ ChannelBlock_ReadAsFrameArrayl32 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

 virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ ChannelBlock_ReadFramesDictUI16 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ ChannelBlock_ReadFramesDict116
 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ ChannelBlock_ReadFramesDictUI32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ ChannelBlock_ReadFramesDictl32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ GetGroupChannelDataUI16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ GetGroupChannelDatal16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ GetGroupChannelDataUI32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< int32_t $>^{\land}>$ $^{\land}$ GetGroupChannelDatal32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

- void SetupGroupDacqQueue (int queuesize, int threshold)
- void SetupGroupDacqQueue (int queuesize, int threshold, unsigned int virtualDevice)
- CHWInfo ^ HWInfo ()

Static Public Attributes

- static const int Error_Callback_Queue_Full = 0x100
- static const int Error Callback Aquisition Stopped = 0x200
- static const int Error_Callback_Packet_Error = 1
- static const int Error_Callback_RingQueue_Full = 3
- static const int Error_Callback_Frames_Lost = 4
- static const int Error_Callback_Data_lost = 5

Properties

virtual int Samplerate [get, set]
 The sampling frequency of the device in Hz.

Events

- OnChannelData^ ChannelDataEvent
- OnError[^] ErrorEvent

Additional Inherited Members

11.48.1 Detailed Description

Base class for data acquisition devices.

11.48.2 Constructor & Destructor Documentation

```
11.48.2.1 CMcsUsbDacqNet() CMcsUsbDacqNet ( )
```

11.48.2.2 ~CMcsUsbDacqNet() ~CMcsUsbDacqNet ()

11.48.3 Member Function Documentation

$\textbf{11.48.3.1} \quad \textbf{AddSelectedChannelsQueue() [1/4]} \quad \text{virtual int AddSelectedChannelsQueue (} \\$

```
int nByteOffset,
int nChannelOffset,
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize ) [virtual]
```

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read ← FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

Parameters

nChannelOffset	Number of channel to start with (counted in samplesize bytes).

Parameters

selectedChannels	List of channels to be collected in the FIFO.
------------------	---

Parameters

queuesize	Size of sample frames the FIFO can hold.
-----------	--

Parameters

threshold	Number of sample frames the FIFO must acquire before the callback function is called.
 แแบบเนา	I MULTIDEL DI SALTIDIE HATTES LIE E IL O HIUSI ACUUITE DEIDTE LIE CAIDACK TUTICIIDITIS CAITCU.

samplesize size of the datawords, either 16 or 32bit.

Returns

The handle to the Queue.

$\textbf{11.48.3.2} \quad \textbf{AddSelectedChannelsQueue() [2/4]} \quad \text{virtual int AddSelectedChannelsQueue (} \\$

```
int nByteOffset,
int nChannelOffset,
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize ) [virtual]
```

11.48.3.3 AddSelectedChannelsQueue() [3/4] virtual int AddSelectedChannelsQueue (

```
int nByteOffset,
int nChannelOffset,
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize ) [virtual]
```

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read ← FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

nByteOffset	Number of bytes to start with.
-------------	--------------------------------

Parameters

Parameters

threshold	Number of sample frames the FIFO must acquire before the callback function is called.

Parameters

```
samplesize size of the datawords, either 16 or 32bit.
```

Returns

The handle to the Queue.

$\textbf{11.48.3.4} \quad \textbf{AddSelectedChannelsQueue() [4/4]} \quad \text{virtual int AddSelectedChannelsQueue (} \\$

```
int nByteOffset,
int nChannelOffset,
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize ) [virtual]
```

```
11.48.3.5 ChannelBlock_AvailFrames() [1/2] virtual uint32_t ChannelBlock_AvailFrames ( int handle ) [virtual]
```

Get the number of sample frames already available in the FIFO.

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.

Returns

Number of sample frames available in the FIFO.

```
11.48.3.6 ChannelBlock_AvailFrames() [2/2] virtual uint32_t ChannelBlock_AvailFrames (
    int handle,
    int queue ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

handle Handle of the FIFO queue. Zero when the SetSelectedData call was used.

Parameters

frames	Number of sample frames to read.
II allies	i Nullibel di Sallible Itallies io lead.

Parameters

frames_ret Number of sample frames which were read, might be smaller than frames.

Returns

Array of int16_t frame arrays.

```
int queue,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Ī

Parameters

queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

Returns

Array of int16_t frame arrays.

Read data from a FIFO queue as array of uint16_t data frame arrays

```
handle Handle of the FIFO queue. Zero when the SetSelectedData call was used.
```

frames Number of sample frames to read
--

Parameters

Returns

Array of int32_t frame arrays.

11.48.3.10 ChannelBlock_ReadAsFrameArrayl32() [2/2] virtual array<array<int32_t>^> ^ Channel← Block_ReadAsFrameArrayI32 (

```
int handle,
int queue,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedData call was used.	
queue	Number of the sub queue.	
frames	Number of sample frames to read.	

Parameters

Returns

Array of int32_t frame arrays.

11.48.3.11 ChannelBlock_ReadAsFrameArrayUl16() [1/2] virtual array<array<uint16_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Parameters

Parameters

frames ret Number of sample frames	which were read, might be smaller than frames.
------------------------------------	--

Returns

Array of uint16_t frame arrays.

11.48.3.12 ChannelBlock_ReadAsFrameArrayUl16() [2/2] virtual array<array<uint16_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

```
handle Handle of the FIFO queue. Zero when the SetSelectedData call was used.
```

queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

Returns

Array of uint16_t frame arrays.


```
ChannelBlock_ReadAsFrameArrayUI32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
--------	--

Parameters

frames	Number of sample frames to read.
--------	----------------------------------

Returns

Array of uint32_t frame arrays.

11.48.3.14 ChannelBlock_ReadAsFrameArrayUl32() [2/2] virtual array<array<uint32_t>^> ^

```
ChannelBlock_ReadAsFrameArrayUI32 (
        int handle,
        int queue,
        int frames,
        [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Parameters

queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

frames_ret	Number of sample frames which were read, might be smaller than frames.
------------	--

Returns

Array of uint32_t frame arrays.

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

dle Handle of the FIFO queue. Zero when the SetSelected	ChannelsQueue call was used.
---	------------------------------

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int16_t arrays and hardware channel as key.

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

handle Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was use
--

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int32_t arrays and hardware channel as key.

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint16_t arrays and hardware channel as key.

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
--------	---

frames Number of sample frames to read.		Number of sample frames to read.
	frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint32_t arrays and hardware channel as key.

```
11.48.3.19 ChannelBlock_ReadFramesI16() [1/2] virtual uint32_t ChannelBlock_ReadFramesI16 (
    int handle,
    array< int16_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int16_t data format

Parameters

handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
--------	---

Parameters

buffer	Buffer to put the data from the device in.	
frames_pos	Position in buffer where to put the data.	
frames	Number of sample frames to read.	

Parameters

frames_ret Number of sample frames which were read, might be smaller than fram
--

Returns

Error Status. 0 on success.

```
11.48.3.20 ChannelBlock_ReadFrames116() [2/2] virtual array<int16_t> ^ ChannelBlock_Read←
FramesI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int16_t data format

	handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.	
--	--------	---	--

Parameters

frames	Number of sample frames to read.
--------	----------------------------------

Parameters

frames_ret Number of sample frames which were read, might be smaller that

```
11.48.3.21 ChannelBlock_ReadFramesI32() [1/2] virtual uint32_t ChannelBlock_ReadFramesI32 (
    int handle,
    array< int32_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.	

Parameters

buffer	Buffer to put the data from the device in.	
frames_pos	Position in buffer where to put the data.	
frames	Number of sample frames to read.	

Returns

Error Status. 0 on success.

```
11.48.3.22 ChannelBlock_ReadFramesl32() [2/2] virtual array<int32_t> ^ ChannelBlock_Read←
FramesI32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.

Parameters

Parameters

frames_ret Number of sample frames which were read, might be smaller than frames.

```
11.48.3.23 ChannelBlock_ReadFramesUl16() [1/2] virtual uint32_t ChannelBlock_ReadFramesUI16 (
    int handle,
    array< uint16_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint16_t data format

handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel	
	number.	

buffer	Buffer to put the data from the device in.
frames_pos	Position in buffer where to put the data.
frames	Number of sample frames to read.

Parameters

frames_ret Number of sample frames v	hich were read, might be smaller than frames.
--------------------------------------	---

Returns

Error Status. 0 on success.

```
11.48.3.24 ChannelBlock_ReadFramesUl16() [2/2] virtual array<uint16_t> ^ ChannelBlock_Read←
FramesUI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint16_t data format

Parameters

handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.
--------	---

Parameters

frames	Number of sample frames to read.
Harries	realiser of sample frames to read.

, ,	
trames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Array of data from the device.

11.48.3.25 ChannelBlock_ReadFramesUI32() [1/2] virtual uint32_t ChannelBlock_ReadFramesUI32 (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out] int % frames_ret) [virtual]

Read data from a FIFO queue in uint32_t data format

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel nun
--

Parameters

buffer	Buffer to put the data from the device in.
frames_pos	Position in buffer where to put the data.
frames	Number of sample frames to read.

Parameters

frames_ret Number of sample frames which were read, might be smaller than fram
--

Returns

Error Status. 0 on success.

```
11.48.3.26 ChannelBlock_ReadFramesUl32() [2/2] virtual array<uint32_t> ^ ChannelBlock_Read←
FramesUl32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the char	iel number.
---	-------------

Parameters

	frames	Number of sample frames to read.
--	--------	----------------------------------

Parameters

frames_ret Number of sample frames which were read, might be smaller than frames.

```
11.48.3.27 ClearBuffers() virtual void ClearBuffers ( ) [virtual]
```

```
11.48.3.28 CMcsUsbDacqNet::GetFilterProperties() virtual array<CFilterPropertyNet^> ^ CMcs↔ UsbDacqNet::GetFilterProperties (

DacqGroupChannelEnumNet GroupID ) [virtual]
```

11.48.3.29 GetAdapterType() virtual AdapterTypeEnumNet GetAdapterType () [virtual]

Gets the adapter which is connected to the MEA2100 device.

Returns

AdapterTypeEnumNet which enumerates the possible adapters.

```
11.48.3.30 GetAdcDataFormat() virtual uint32_t GetAdcDataFormat ( uint32_t virtualDevice ) [virtual]
```

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

Returns

The data format in bits.

```
11.48.3.31 GetAdcZero() virtual uint32_t GetAdcZero (
             uint32_t virtualDevice,
             DacqGroupChannelEnumNet group,
             [System::Runtime::InteropServices::Out] int% adcz ) [virtual]
11.48.3.32 GetAnalogValueUnit() virtual uint32_t GetAnalogValueUnit (
             uint32_t virtualDevice,
             DacqGroupChannelEnumNet group,
             [System::Runtime::InteropServices::Out] AnalogUnitEnumNet% unit ) [virtual]
11.48.3.33 GetChannelDataFillSize() int GetChannelDataFillSize ()
11.48.3.34 GetChannelLayout() virtual uint32_t GetChannelLayout (
             [System::Runtime::InteropServices::Out] int% AnalogChannels,
             [System::Runtime::InteropServices::Out] int% DigitalChannels,
             [System::Runtime::InteropServices::Out] int% ChecksumChannels,
             [System::Runtime::InteropServices::Out] int% TimestampChannels,
             [{\tt System::Runtime::InteropServices::Out}] \  \, {\tt int} \$ \  \, {\tt ChannelsInBlock},
             unsigned int virtualDevice ) [virtual]
11.48.3.35 GetChannelsInBlock() virtual int32_t GetChannelsInBlock (
             unsigned int virtualDevice ) [virtual]
```

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

Returns

Number of 16 bit datawords per sample frame.

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

virtualDevice	Virtual device to use.
---------------	------------------------

Returns

DataModeEnumNet which enumerates the possible data modes.

Gets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the MEA2100 device.

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber offset	An offset / bit number with the source/function.

[System::Runtime::InteropServices::Out] int% bitnumber_offset)

Gets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the W2100 device.

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber offset Generated by Doxygen	An offset / bit number with the source/function.

```
11.48.3.43 GetFilterProperty() virtual CFilterPropertyNet ^ GetFilterProperty (
DacqGroupChannelEnumNet GroupID,
unsigned int index ) [virtual]
```

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

elector supported by the device.	group Group selector
----------------------------------	----------------------

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int16_t arrays and hardware channel as key.

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

ı		
	group	Group selector supported by the device.

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int32_t arrays and hardware channel as key.

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

	Outside a state of a state of the state of t
group	Group selector supported by the device.

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint16_t arrays and hardware channel as key.

```
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

group	Group selector supported by the device.
-------	---

Parameters

	frames	Number of sample frames to read.
	frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint32 t arrays and hardware channel as key.

```
11.48.3.48 GetHardwareMaxRange() virtual uint32_t GetHardwareMaxRange ( uint32_t virtualDevice,
```

```
DacqGroupChannelEnumNet group,
[System::Runtime::InteropServices::Out] int% r,
[System::Runtime::InteropServices::Out] int% rUnit ) [virtual]
```

11.48.3.49 GetHardwareMinRange() virtual uint32_t GetHardwareMinRange (

```
uint32_t virtualDevice,
DacqGroupChannelEnumNet group,
[System::Runtime::InteropServices::Out] int% r,
[System::Runtime::InteropServices::Out] int% rUnit ) [virtual]
```

```
11.48.3.50 GetMaxSamplingFrequency() virtual uint32_t GetMaxSamplingFrequency ( int virtualDevice ) [virtual]
```

Gets the maximal sampling frequency of the device.

Returns

Sampling frequency in Hz.

```
11.48.3.51 GetMeaLayout() virtual MeaLayoutEnumNet GetMeaLayout ( ) [virtual]
```

Gets the MEA layout which is connected to the MEA2100 device.

Returns

MeaLayoutEnumNet which enumerates the MEA types.

```
11.48.3.52 GetMinSamplingFrequencyStepsize() virtual uint32_t GetMinSamplingFrequencyStepsize
( ) [virtual]
```

Gets the minimal sampling frequency step size increment value of the device.

Returns

Sampling frequency step size in Hz.

Get the real number of data bits.

This value may be different from the value returned by GetDataFormat, e.g. in MC_Card the data are shifted 2 bits so the real number is 14 while the data format is 16 bits

```
11.48.3.56 GetSamplerate() virtual int32_t GetSamplerate ( unsigned int virtualDevice ) [virtual]
```

Gets the sampling frequency of the device.

Returns

Sampling frequency in Hz.

```
\textbf{11.48.3.57} \quad \textbf{GetVoltageRangeIndex()} \quad \texttt{virtual uint} \\ \textbf{12\_t GetVoltageRangeIndex} \ (
              unsigned int virtualDevice ) [virtual]
11.48.3.58 GetVoltageRangeInMicroVolt() virtual int32_t GetVoltageRangeInMicroVolt (
              unsigned int virtualDevice ) [virtual]
Gets the currently selected voltage range on devices which support multiple voltage ranges.
Returns
     The Voltage Range in uV.
11.48.3.59 GetVoltageRangeInMilliVolt() virtual int32_t GetVoltageRangeInMilliVolt ( ) [virtual]
Gets the currently selected voltage range on devices which support multiple voltage ranges.
Returns
     The rounded Voltage Range in mV.
11.48.3.60 HWInfo() CHWInfo ^ HWInfo ( )
11.48.3.61 SendStartDacq() [1/2] virtual void SendStartDacq ( ) [virtual]
Start sampling.
11.48.3.62 SendStartDacq() [2/2] virtual void SendStartDacq (
              int VirtualDacqMap ) [virtual]
Start sampling.
Parameters
 VirtualDacqMap
```

Start sampling together with the STG.

Parameters

trigger_map

VirtualDacqMap

11.48.3.64 SendStopDacq() [1/2] virtual void SendStopDacq () [virtual]

Stop sampling.

Stop sampling.

Parameters

VirtualDacqMap

Stop sampling together with the STG.

Parameters

trigger_map

```
11.48.3.67 SendStopStgAndDacqWithOptions() virtual void SendStopStgAndDacqWithOptions ( uint32_t trigger_map,
```

```
int VirtualDacqMap,
int options ) [virtual]
```

Stop sampling together with the STG and options.

Parameters

trigger_map

Parameters

options

Parameters

VirtualDacqMap

```
11.48.3.68 SetDataMode() virtual void SetDataMode (
DataModeEnumNet dataMode,
unsigned int virtualDevice ) [virtual]
```

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

dataMode	DataModeEnumNet enumerates the possible data modes.
virtualDevice	Virtual device to use.

Sets the function/source of an digital output bit.

This is the templated generic implementation.

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

```
11.48.3.72 SetDigitalSource() [4/5] void SetDigitalSource (
DigitalTargetEnumNet digitaltarget,
```

```
int32_t NrChannel,
SCUDigitalSourceEnumNet source,
int bitnumber_offset )
```

Sets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the sampling frequency of the device.

rate Sampling frequency in Hz.

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be devided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

selectedChannels	List of channels to be collected in the FIFO.
------------------	---

Parameters

queuesize	Size of sample frames the FIFO can hold.
-----------	--

Parameters

t	hreshold	Number of sample	frames the FIFO must	acquire before the	callback function is called.
---	----------	------------------	----------------------	--------------------	------------------------------

Parameters

samplesize	size of the datawords, either 16 or 32bit.
------------	--

11.48.3.77 SetSelectedChannels() [2/4] virtual void SetSelectedChannels (

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

11.48.3.78 SetSelectedChannels() [3/4] virtual void SetSelectedChannels (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock ) [virtual]
```

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be devided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

nChannele	Number of channels to be collected in the FIFO.
HUHAHHUS	i Nullibel di chalifiels lo de collected ili tile i il O.

Parameters

queuesize S	Size of sample frames the FIFO can hold.
-------------	--

Parameters

throchold	Number of samples frames the FIFO must acquire before the callback function is called.
unesnou	Number of samples frames the Firo must acquire before the caliback function is called.

samplesize	size of the datawords, either 16 or 32bit.				
ChannelsInBlock	value obtained from GetChannelsInBlock.				

11.48.3.79 SetSelectedChannels() [4/4] virtual void SetSelectedChannels (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

11.48.3.80 SetSelectedChannelsQueue() [1/4] virtual void SetSelectedChannelsQueue (

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock ) [virtual]
```

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

selectedChannels	List of channels to be collected in the FIFO.

Parameters

queuesize Size of sample frames the FIFO	can hold.
--	-----------

Parameters

threshold	Number of sample frames the FIFO must acquire before the callback function is called.
-----------	---

	aine af the determent aither 10 an 00bit
sambiesize	size of the datawords, either 16 or 32bit.

```
11.48.3.81 SetSelectedChannelsQueue() [2/4] virtual void SetSelectedChannelsQueue (
```

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock) [virtual]
```

11.48.3.82 SetSelectedChannelsQueue() [3/4] virtual void SetSelectedChannelsQueue (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock ) [virtual]
```

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

nChannels	Number of channels to be collected in the FIFO.
-----------	---

Parameters

queuesize	Size of sample frames the FIFO can hold.
-----------	--

thrac	hald I	Numbero	f sample fram	oc tha FIF∩	must acquire	hafora tha	callback fui	nction ic	called !
unco	ioia	I NUITIDEI O	i sample mam		illust acquire		caliback ful		canca.

samplesize size of the datawords, either 16 or 32bit.

Parameters

ChannelsInBlock	value obtained from GetChannelsInBlock.
Unanifelaniblock	value obtained nom detonameismblock

11.48.3.83 SetSelectedChannelsQueue() [4/4] virtual void SetSelectedChannelsQueue (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

11.48.3.84 SetSelectedData() [1/4] virtual void SetSelectedData (

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock) [virtual]
```

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

selectedChannels	List of channels to be collected in the FIFO.
------------------	---

auguggiza	Size of sample frames the FIFO can hold.
uueuesize	. Size of Samble frames the FIFO Can hold.

threshold	Number of sample frames the FIFO must acquire before the callback function is called.
-----------	---

Parameters

samplesize	size of the datawords, either 16 or 32bit.
ChannelsInBlock	value obtained from GetChannelsInBlock.

int ChannelsInBlock) [virtual]

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

nChannels	Number of channels to be collected in the FIFO.

queuesize	Size of sample frames the FIFO can hold.

Parameters

samplesize	size of the datawords, either 16 or 32bit.	
ChannelsInBlock	value obtained from GetChannelsInBlock.	

```
11.48.3.87 SetSelectedData() [4/4] virtual void SetSelectedData (
    int nChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    SampleDstSizeNet sampleDstSize,
    int ChannelsInBlock ) [virtual]

11.48.3.88 SetupGroupDacqQueue() [1/2] void SetupGroupDacqQueue (
    int queuesize,
    int threshold)

11.48.3.89 SetupGroupDacqQueue() [2/2] void SetupGroupDacqQueue (
    int queuesize,
    int threshold,
    unsigned int virtualDevice )

11.48.3.90 SetVoltageRangeByIndex() virtual void SetVoltageRangeByIndex (
```

Sets the voltage range on devices which support multiple voltage ranges.

unsigned int virtualDevice) [virtual]

int32_t voltageRangeIndex,

		- 1
	Voltage Range to use as index, smaller values are larger voltage ranges.	- 1
VoltadeKandelndey I	. Voltage Bange to like as index, smaller vallies are larger voltage ranges	- 1
voitager larigerilaex	i voltage i lange to use as index, sindher values are langer voltage ranges.	- 1

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

voltageRange	Voltage Range to use in μV.
1 0 110.0 9 0 1 10.11 9 0	i consign consign to see in pro-

This replaces SetVoltageRange, where the value of the range was in mV!

```
11.48.3.92 StartDacq() [1/4] virtual void StartDacq ( ) [virtual]
```

Start the data acquisition thread and sampling.

```
11.48.3.93 StartDacq() [2/4] virtual void StartDacq ( int32_t timeout ) [virtual]
```

Start the data acquisition thread and sampling.

Parameters

```
timeout Timeout in ms.
```

Start the data acquisition thread and sampling.

timeout Time	eout in ms.
--------------	-------------

mber of USB Buffers that are simultaniousl	numSubmittedUsbBuffers
--	------------------------

Parameters

Parameters

packetsInUrb	Packets in each URB.
--------------	----------------------

```
11.48.3.95 StartDacq() [4/4] virtual void StartDacq (
```

```
int32_t timeout,
int32_t numSubmittedUsbBuffers,
int32_t numUsbBuffers,
int32_t packetsInUrb,
uint32_t virtualDevice) [virtual]
```

Start the data acquisition thread and sampling.

Parameters

numSubmittedUsbBuffers	Number of USB Buffers that are simultaniously submitted
------------------------	---

Parameters

timeout	Timeout in ms.
unicuu	Tillicout III IIIo.

Parameters

virtualDevice Vir	tual Device to start.
-------------------	-----------------------

11.48.3.96 StartLoop() [1/4] virtual void StartLoop () [virtual]

Start the data acquisition thread.

Start the data acquisition thread.

Parameters

```
timeout | Timeout in ms.
```

int32_t packetsInUrb) [virtual]

_ 1

Start the data acquisition thread.

Parameters

timeout	Timeout in ms.

numSubmittedUsbBuffers Number of USB Buffers that are simultaniously submitted
--

numUsbBuffers Number of USB Buffe	ers to use.
-----------------------------------	-------------

Parameters

packetsInUrb	Packets in each URB.
--------------	----------------------

```
11.48.3.99 StartLoop() [4/4] virtual void StartLoop() int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_t packetsInUrb, uint32_t virtualDevice) [virtual]
```

Start the data acquisition thread.

Parameters

numSubmittedLlehBuffere	Number of USB Buffers that are simultaniously submitted.
I HUHISUDHIILLEUUSDDUHEIS	I Mullipel of OSD bullers that are simulationary submitted.

Parameters

timeout	Timeout in ms.
---------	----------------

Parameters

	1
millah Duffara	Number of USB Buffers to use.
numusnauners	Number of USB Bullers to USE

packetsInUrb	Packets in each URB.
pachetallield	i donoto ili cdoli olib.

virtualDevice Virtual Device to start.
--

11.48.3.100 StopDacq() [1/2] virtual void StopDacq () [virtual]

Stop the data acquisition thread and sampling.

```
11.48.3.101 StopDacq() [2/2] virtual void StopDacq ( uint32_t virtualDevice ) [virtual]
```

Stop the data acquisition thread and sampling.

Parameters

virtualDevice Virtual Device to start

11.48.3.102 StopLoop() virtual void StopLoop () [virtual]

11.48.4 Member Data Documentation

11.48.4.1 Error_Callback_Aquisition_Stopped const int Error_Callback_Aquisition_Stopped = 0x200 [static]

11.48.4.2 Error_Callback_Data_lost const int Error_Callback_Data_lost = 5 [static]

11.48.4.3 Error_Callback_Frames_Lost const int Error_Callback_Frames_Lost = 4 [static]

11.48.4.4 Error_Callback_Packet_Error const int Error_Callback_Packet_Error = 1 [static]

11.48.4.5 Error_Callback_Queue_Full const int Error_Callback_Queue_Full = 0x100 [static]

11.48.4.6 Error_Callback_RingQueue_Full const int Error_Callback_RingQueue_Full = 3 [static]

11.48.5 Property Documentation

11.48.5.1 Samplerate virtual int Samplerate [get], [set]

The sampling frequency of the device in Hz.


11.48.6 Event Documentation

11.48.6.1 ChannelDataEvent OnChannelData^ ChannelDataEvent

11.48.6.2 ErrorEvent OnError^ ErrorEvent

11.49 CMcsUsbDeviceStatePushFunctionNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushFunctionNet:

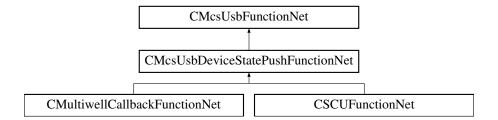

Public Member Functions

• void TriggerStatus ()

Protected Member Functions

CMcsUsbDeviceStatePushFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] p

 Device)

Events

• OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add, remove, raise]

Additional Inherited Members

11.49.1 Constructor & Destructor Documentation

```
11.49.1.1 CMcsUsbDeviceStatePushFunctionNet() CMcsUsbDeviceStatePushFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pDevice ) [protected]
```

11.49.2 Member Function Documentation


```
11.49.2.1 TriggerStatus() void TriggerStatus ()
```

11.49.3 Event Documentation

```
11.49.3.1 McsUsbDeviceStateEvent OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]
```

11.50 CMcsUsbDeviceStatePushNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushNet:

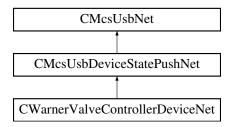

Public Member Functions

• void TriggerStatus ()

Protected Member Functions

CMcsUsbDeviceStatePushNet (CMcsUsbPointerContainer[∧] pDevice)

Events

• OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add, remove, raise]

Additional Inherited Members

11.50.1 Constructor & Destructor Documentation

```
11.50.1.1 CMcsUsbDeviceStatePushNet() CMcsUsbDeviceStatePushNet (
CMcsUsbPointerContainer^ pDevice) [protected]
```

11.50.2 Member Function Documentation


```
11.50.2.1 TriggerStatus() void TriggerStatus ()
```

11.50.3 Event Documentation

```
11.50.3.1 McsUsbDeviceStateEvent OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]
```

11.51 CMcsUsbFactoryNet Class Reference

Inheritance diagram for CMcsUsbFactoryNet:

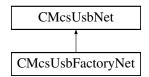

Public Member Functions

- CMcsUsbFactoryNet ()
- ∼CMcsUsbFactoryNet ()
- unsigned int GetNumDestinations ()
- String \(^{\text{GetDestinationName}}\) (unsigned int index)
- String \(^\) GetDestinationName (CFirmwareDestinationNet dest)
- void SetDestinationSerialNumber (CFirmwareDestinationNet dest, String[^] serialnumber)
- String \(^\) GetDestinationSerialNumber (CFirmwareDestinationNet dest)
- CFirmwareDestinationNet GetDestination (unsigned int index)
- CFirmwareDestinationNet GetDestination (String[^] Key)
- unsigned int GetDestinationTargetAddress (CFirmwareDestinationNet destination)

Gets the target base address for the destination.

- uint32_t ChangeSerialNumber (String^ serial)
- bool LoadUserFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry)

Send the DSP Firmware to the MEA21 device.

- bool LoadUserFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, uint32 t LockMask)
- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, bool SkipWait)

Flashes a firmware file to the device.

- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, bool SkipWait, unsigned int LockMask)
- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest)

Flashes a firmware file to the device.

 bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest, bool SkipWait)

Flashes a firmware file to the device.

- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest, bool SkipWait, uint32_t LockMask)
- bool CompareFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestination
 — Net Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, String[^]
 MessagePrefix, unsigned int LockMask, [System::Runtime::InteropServices::Out] String[^]% ErrorText,
 [System::Runtime::InteropServices::Out] String[^]% Protokoll)
- uint32_t Coldstart (CFirmwareDestinationNet dest)
- int32 t GetXilinxFlashOffset (CFirmwareDestinationNet dest)
- uint32 t GetXilinxFlashReadCommand (CFirmwareDestinationNet dest)
- array< uint8 t > ^ DownloadFirmware (CFirmwareDestinationNet Dest, uint32 t Address, uint32 t length)
- bool GetUsercodeFromFlash (unsigned int FPGA, unsigned int Address, [System::Runtime::Interop
 — Services::Out] unsigned int% Usercode)
- array< unsigned char > ^ ReadBlockFromFlash (unsigned int FPGA, unsigned int Address)
- void ReadBlockFromFlash (unsigned int FPGA, unsigned int Address, array< unsigned char >^ buffer, int position)
- array< unsigned char > ^ ReadBlockFromIFBGlobalEEprom (unsigned int Address)
- array< unsigned char > ^ ReadBlockFromNVMEM (unsigned int FPGA, unsigned int Offset, unsigned int Address)

Static Public Member Functions

- static String \(^\) GetDestinationDisplayLabel (String\(^\) RawLabel, CFirmwareDestinationNet dest)
- static String ^ FindFirmwareVersionMagicInBuffer (array< unsigned char >^ buffer, int length, [System::

 Runtime::InteropServices::Out]int% position)
- static bool GetFirmwareVersionFromFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] uint32 t% Version)

Retrives version info from a Firmware update file.

- static bool GetFirmwareVersionFromFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version, [System::Runtime::InteropServices::Out] uint32_t% Position)
- static bool GetFirmwareVersionFromHexFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
- static uint32_t GetChecksumFromFX3Image (String[^] FirmwareFile)
- static uint32_t GetUSBDeviceIDFromFX3Image (String[^] FirmwareFile)
- static bool GetUsercodeFromBitFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] unsigned int% Usercode)

Static Public Attributes

- static const uint32 t FX3MCSDataAddress = 0x40037E00
- static const uint32_t FX3MCSDataDeviceIdOffset = 0x4
- static const uint32_t FX3MCSDataVersionOffset = 0x8
- static const uint32_t FX3MCSDataIFB2ImageOffset = 0xC
- static const uint32_t FX3MCSDataIFB1ImageOffset = 0x2C

Additional Inherited Members

11.51.1 Constructor & Destructor Documentation

```
11.51.1.1 CMcsUsbFactoryNet() CMcsUsbFactoryNet ( )
```

11.51.1.2 ~CMcsUsbFactoryNet() ~CMcsUsbFactoryNet ()

11.51.2 Member Function Documentation

```
11.51.2.1 ChangeSerialNumber() uint32_t ChangeSerialNumber ( String^{\wedge} serial)
```

```
11.51.2.2 Coldstart() uint32_t Coldstart (
              CFirmwareDestinationNet dest)
11.51.2.3 CompareFirmware() bool CompareFirmware (
              String^ FirmwareFile,
              CMcsUsbListEntryNet^ listEntry,
              CFirmwareDestinationNet Dest,
              OnUpdateFirmwareStatusChange^ deleg,
              OnUpdateFirmwareProgress^ progress,
              String^{\wedge} MessagePrefix,
              unsigned int LockMask,
              [System::Runtime::InteropServices::Out] String^% ErrorText,
              [System::Runtime::InteropServices::Out] String^% Protokoll )
11.51.2.4 DownloadFirmware() array<uint8_t> ^ DownloadFirmware (
              CFirmwareDestinationNet Dest,
              uint32_t Address,
              uint32_t length )
11.51.2.5 FindFirmwareVersionMagicInBuffer() static String ^ FindFirmwareVersionMagicInBuffer (
              array< unsigned char >^{\wedge} buffer,
              int length,
              [System::Runtime::InteropServices::Out] int% position ) [static]
11.51.2.6 GetChecksumFromFX3Image() static uint32_t GetChecksumFromFX3Image (
              String<sup>∧</sup> FirmwareFile ) [static]
11.51.2.7 GetDestination() [1/2] CFirmwareDestinationNet GetDestination (
              String^{\wedge} Key )
11.51.2.8 GetDestination() [2/2] CFirmwareDestinationNet GetDestination (
              unsigned int index)
\textbf{11.51.2.9} \quad \textbf{GetDestinationDisplayLabel()} \quad \texttt{static String} \; \land \; \texttt{GetDestinationDisplayLabel} \; \; \textbf{(}
              String^ RawLabel,
              CFirmwareDestinationNet dest ) [static]
```

```
11.51.2.10 GetDestinationName() [1/2] String ^ GetDestinationName (
             CFirmwareDestinationNet dest )
11.51.2.11 GetDestinationName() [2/2] String ^ GetDestinationName (
             unsigned int index )
11.51.2.12 GetDestinationSerialNumber() String ^ GetDestinationSerialNumber (
             CFirmwareDestinationNet dest)
11.51.2.13 GetDestinationTargetAddress() unsigned int GetDestinationTargetAddress (
             CFirmwareDestinationNet destination )
Gets the target base address for the destination.
Parameters
 destination
             The destination to be queried.
Returns
     The base address as a 32 bit number, only the lower 16 bit represent the address.
11.51.2.14 GetFirmwareVersionFromFile() [1/2] static bool GetFirmwareVersionFromFile (
             String^ FirmwareFile,
             [System::Runtime::InteropServices::Out] uint32_t% Version ) [static]
Retrives version info from a Firmware update file.
11.51.2.15 GetFirmwareVersionFromFile() [2/2] static bool GetFirmwareVersionFromFile (
             String Firmware File,
             [System::Runtime::InteropServices::Out] uint32_t% Version,
             [System::Runtime::InteropServices::Out] uint32_t% Position ) [static]
11.51.2.16 GetFirmwareVersionFromHexFile() static bool GetFirmwareVersionFromHexFile (
```

[System::Runtime::InteropServices::Out] uint32_t% Version) [static]

String^ FirmwareFile,

```
11.51.2.17 GetNumDestinations() unsigned int GetNumDestinations ( )
\textbf{11.51.2.18} \quad \textbf{GetUSBDeviceIDFromFX3Image()} \quad \texttt{static uint32\_t GetUSBDeviceIDFromFX3Image ()} \\
               String FirmwareFile ) [static]
11.51.2.19 GetUsercodeFromBitFile() static bool GetUsercodeFromBitFile (
               String<sup>^</sup> FirmwareFile,
               [{\tt System::Runtime::InteropServices::Out}] \ \ unsigned \ \ int \ \ \textit{Usercode} \ ) \ \ \ [{\tt static}]
\textbf{11.51.2.20} \quad \textbf{GetUsercodeFromFlash()} \quad \texttt{bool GetUsercodeFromFlash ()}
               unsigned int FPGA,
               unsigned int Address,
               [System::Runtime::InteropServices::Out] \ unsigned \ int \% \ \textit{Usercode} \ )
11.51.2.21 GetXilinxFlashOffset() int32_t GetXilinxFlashOffset (
               CFirmwareDestinationNet dest )
11.51.2.22 GetXilinxFlashReadCommand() uint32_t GetXilinxFlashReadCommand (
               CFirmwareDestinationNet dest )
11.51.2.23 LoadUserFirmware() [1/2] bool LoadUserFirmware (
               String^{\wedge} FirmwareFile,
               CMcsUsbListEntryNet^ listEntry )
Send the DSP Firmware to the MEA21 device.
Parameters
 FirmwareFile
                 Filename of the DSP Firmware (*.bin) file.
```

listEntry Device to use for the connection. See CMcsUsbListNet.

```
11.51.2.24 LoadUserFirmware() [2/2] bool LoadUserFirmware (
             String^ FirmwareFile,
             CMcsUsbListEntryNet^ listEntry,
             uint32_t LockMask )
11.51.2.25 ReadBlockFromFlash() [1/2] array<unsigned char> ^ ReadBlockFromFlash (
             unsigned int FPGA,
             unsigned int Address )
11.51.2.26 ReadBlockFromFlash() [2/2] void ReadBlockFromFlash (
             unsigned int FPGA,
             unsigned int Address,
             array< unsigned char >^{\wedge} buffer,
             int position )
11.51.2.27 ReadBlockFromIFBGlobalEEprom() array<unsigned char> ^ ReadBlockFromIFBGlobalE←
Eprom (
             unsigned int Address )
11.51.2.28 ReadBlockFromNVMEM() array<unsigned char> ^ ReadBlockFromNVMEM (
             unsigned int FPGA,
             unsigned int Offset,
             unsigned int Address )
11.51.2.29 SetDestinationSerialNumber() void SetDestinationSerialNumber (
             CFirmwareDestinationNet dest,
             String^ serialnumber )
11.51.2.30 UpdateFirmware() [1/5] bool UpdateFirmware (
             String^ FirmwareFile,
             CMcsUsbListEntryNet<sup>^</sup> listEntry,
             CFirmwareDestinationNet dest )
```

Flashes a firmware file to the device.

FirmwareFile	Filename of the Firmware file.	
--------------	--------------------------------	--

Parameters

listEntry Device to use for the co	onnection.
------------------------------------	------------

Flashes a firmware file to the device.

Parameters

Parameters

```
listEntry Device to use for the connection.
```

bool SkipWait,
uint32_t LockMask)

```
\textbf{11.51.2.33} \quad \textbf{UpdateFirmware() [4/5]} \quad \texttt{bool UpdateFirmware (}
```

```
String FirmwareFile,
CMcsUsbListEntryNet listEntry,
CFirmwareDestinationNet Dest,
```

```
OnUpdateFirmwareStatusChange^ deleg,
OnUpdateFirmwareProgress^ progress,
bool SkipWait )
```

Flashes a firmware file to the device.

Parameters

```
FirmwareFile Filename of the Firmware file.
```

```
11.51.2.34 UpdateFirmware() [5/5] bool UpdateFirmware (
String^ FirmwareFile,
```

```
CMcsUsbListEntryNet<sup>^</sup> listEntry,
CFirmwareDestinationNet Dest,
OnUpdateFirmwareStatusChange<sup>^</sup> deleg,
OnUpdateFirmwareProgress<sup>^</sup> progress,
bool SkipWait,
unsigned int LockMask)
```


11.51.3 Member Data Documentation

```
11.51.3.1 FX3MCSDataAddress const uint32_t FX3MCSDataAddress = 0x40037E00 [static]
```

- 11.51.3.2 FX3MCSDataDeviceIdOffset const uint32_t FX3MCSDataDeviceIdOffset = 0x4 [static]
- 11.51.3.3 FX3MCSDataIFB1ImageOffset const uint32_t FX3MCSDataIFB1ImageOffset = 0x2C [static]
- 11.51.3.4 FX3MCSDataIFB2ImageOffset const uint32_t FX3MCSDataIFB2ImageOffset = 0xC [static]
- 11.51.3.5 FX3MCSDataVersionOffset const uint32_t FX3MCSDataVersionOffset = 0x8 [static]

11.52 CMcsUsbFunctionNet Class Reference

Inheritance diagram for CMcsUsbFunctionNet:

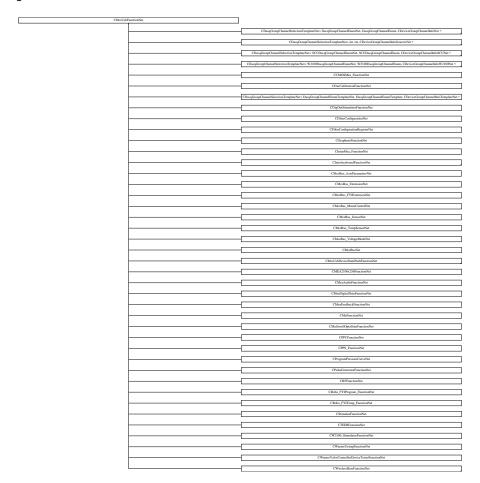

Public Member Functions

- CMcsUsbFunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual \sim CMcsUsbFunctionNet (void)
- !CMcsUsbFunctionNet ()
- void ThrowCUsbExceptionNetOnError (uint32_t status)

Protected Member Functions

• CMcsUsbFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] mcsusbfunction)

Protected Attributes

- CMcsUsbNet ^ m_pMcsUsb
- CMcsUsbFunction * m_pMcsUsbFunction

11.52.1 Constructor & Destructor Documentation

```
11.52.1.1 CMcsUsbFunctionNet() [1/2] CMcsUsbFunctionNet (
            CMcsUsbNet^ mcsusb )
11.52.1.2 ~CMcsUsbFunctionNet() virtual ~CMcsUsbFunctionNet (
            void ) [virtual]
11.52.1.3 "!CMcsUsbFunctionNet() !CMcsUsbFunctionNet ()
11.52.1.4 CMcsUsbFunctionNet() [2/2] CMcsUsbFunctionNet (
            CMcsUsbNet^ mcsusb,
            CMcsUsbFunctionPointerContainer<sup>∧</sup> mcsusbfunction ) [protected]
11.52.2 Member Function Documentation
11.52.2.1 ThrowCUsbExceptionNetOnError() void ThrowCUsbExceptionNetOnError (
            uint32_t status )
11.52.3 Member Data Documentation
11.52.3.1 m pMcsUsb CMcsUsbNet ^ m_pMcsUsb [protected]
11.52.3.2 m pMcsUsbFunction CMcsUsbFunction* m_pMcsUsbFunction [protected]
11.53 CMcsUsbFunctionPointerContainer Class Reference
11.54 CMcsUsbListEntryNet Class Reference
```

Generated by Doxygen

McsUsbListEntryNet identifies a connected device.

Public Member Functions

- ∼CMcsUsbListEntryNet ()
- virtual bool Equals (Object[^] obj) override

Checks weather two CMcsUsbListEntryNet represent the same USB device.

void SetStringFormat (String ^ format)

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

• virtual String ^ ToString () override

Static Public Member Functions

static CMcsUsbListEntryNet [^] GetEntry ()

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

static CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

static CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice, unsigned int index)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

static unsigned int GetEntryCount ()

Returns the number of devices connected to the computer.

static unsigned int GetEntryCount (DeviceEnumNet McsUsbDevice)

Returns the number of devices connected to the computer.

Properties

• String Manufacturer [get]

The Manufacturer ID of the device represented by this CMcsUsbListEntryNet.

String[^] Product [get]

The Product ID of the device represented by this CMcsUsbListEntryNet.

• String DeviceName [get]

The device name of the device represented by this CMcsUsbListEntryNet.

String[^] SerialNumber [get]

The serial number of the device represented by this CMcsUsbListEntryNet.

• String HwVersion [get]

The hardware revision of the device represented by this CMcsUsbListEntryNet.

• DeviceIdNet^ DeviceId [get]

11.54.1 Detailed Description

McsUsbListEntryNet identifies a connected device.

11.54.2 Constructor & Destructor Documentation

11.54.2.1 ~CMcsUsbListEntryNet() ~CMcsUsbListEntryNet ()

11.54.3 Member Function Documentation

```
11.54.3.1 Equals() virtual bool Equals (
Object^ obj ) [override], [virtual]
```

Checks weather two CMcsUsbListEntryNet represent the same USB device.

obj The CMcsUsbListEntryNet to compare with.

```
11.54.3.2 GetEntry() [1/3] static CMcsUsbListEntryNet ^ GetEntry ( ) [static]
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

```
11.54.3.3 GetEntry() [2/3] static CMcsUsbListEntryNet ^ GetEntry (

DeviceEnumNet McsUsbDevice ) [static]
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

McsUsbDevice | Specifies the type of devices to look for.

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

McsUsbDevice | Specifies the type of devices to look for.

	index	number of the entry to use.
--	-------	-----------------------------

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

11.54.3.5 GetEntryCount() [1/2] static unsigned int GetEntryCount () [static]

Returns the number of devices connected to the computer.

Returns

The number of devices.

```
11.54.3.6 GetEntryCount() [2/2] static unsigned int GetEntryCount (

DeviceEnumNet McsUsbDevice ) [static]
```

Returns the number of devices connected to the computer.

Parameters

Returns

The number of devices.

```
11.54.3.7 SetStringFormat() void SetStringFormat (
String ^ format )
```

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

format A String containing the format template. Default is "%N (%S)".

11.54.3.8 ToString() virtual String ^ ToString () [override], [virtual]

11.54.4 Property Documentation

11.54.4.1 DeviceId DeviceIdNet^ DeviceId [get]

11.54.4.2 DeviceName String DeviceName [get]

The device name of the device represented by this CMcsUsbListEntryNet.

11.54.4.3 HwVersion String MwVersion [get]

The hardware revision of the device represented by this CMcsUsbListEntryNet.

11.54.4.4 Manufacturer String^ Manufacturer [get]

The Manufacturer ID of the device represented by this CMcsUsbListEntryNet.

11.54.4.5 **Product** String^ Product [get]

The Product ID of the device represented by this CMcsUsbListEntryNet.

11.54.4.6 SerialNumber String SerialNumber [get]

The serial number of the device represented by this CMcsUsbListEntryNet.

11.55 CMcsUsbListNet Class Reference

Class to handle a list of connected MCS USB devices.

Public Member Functions

CMcsUsbListNet (DeviceEnumNet McsUsbDevice)

Initializes a new instance of CMcsUsbListNet class.

CMcsUsbListNet (array< DeviceIdNet[^]>[^] DeviceIdList)

Initializes a new instance of CMcsUsbListNet class.

∼CMcsUsbListNet ()

Destructor: called by Dispose()

• !CMcsUsbListNet ()

Finalizer: called by GC before collecting

void SetStringFormat (String ^ format)

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

• uint32 t GetNumberOfDevices ()

Gets the number of devices currently in the list.

CMcsUsbListEntryNet ^ GetUsbListEntry (unsigned int index)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

array< CMcsUsbListEntryNet[^]> [^] GetUsbListEntries ()

Returns all entries from the list of USB Devices connected to the computer.

bool IsDeviceTypeOf (CMcsUsbListEntryNet[^] entry, DeviceEnumNet McsUsbDevice)

Properties

• uint32_t Count [get]

Gets the number of devices currently in the list.

Events

- OnDeviceArrivalRemoval[^] DeviceArrival
- OnDeviceArrivalRemoval[^] DeviceRemoval

11.55.1 Detailed Description

Class to handle a list of connected MCS USB devices.

11.55.2 Constructor & Destructor Documentation

```
11.55.2.1 CMcsUsbListNet() [1/2] CMcsUsbListNet (

DeviceEnumNet McsUsbDevice)
```

Initializes a new instance of CMcsUsbListNet class.

Initializes a new instance of CMcsUsbListNet class.

```
11.55.2.3 ~CMcsUsbListNet() ~CMcsUsbListNet ()
```

Destructor: called by Dispose()

11.55.2.4 "!CMcsUsbListNet() !CMcsUsbListNet ()

Finalizer: called by GC before collecting

11.55.3 Member Function Documentation

```
11.55.3.1 GetNumberOfDevices() uint32_t GetNumberOfDevices ( )
```

Gets the number of devices currently in the list.

Returns

The number of devices currently in the list.

```
11.55.3.2 GetUsbListEntries() array<CMcsUsbListEntryNet^> ^ GetUsbListEntries ( )
```

Returns all entries from the list of USB Devices connected to the computer.

```
11.55.3.3 GetUsbListEntry() CMcsUsbListEntryNet ^ GetUsbListEntry ( unsigned int index )
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

index | number of the entry to use.

```
11.55.3.5 SetStringFormat() void SetStringFormat ( String ^{\wedge} format )
```

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

format A String containing the format template. Default is "%	N (%S)".
---	----------

11.55.4 Property Documentation

11.55.4.1 Count uint32_t Count [get]

Gets the number of devices currently in the list.

11.55.5 Event Documentation


11.55.5.1 DeviceArrival OnDeviceArrivalRemoval^ DeviceArrival

11.55.5.2 DeviceRemoval OnDeviceArrivalRemoval^ DeviceRemoval

11.56 CMcsUsbNet Class Reference

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

Inheritance diagram for CMcsUsbNet:

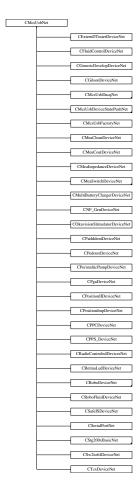

Public Member Functions

• CMcsUsbNet ()

Initializes a new instance of the base class to handle MCS USB devices.

CMcsUsbNet (McsBusTypeEnumNet bustype)

Initializes a new instance of the base class to handle MCS USB devices.

- virtual ∼CMcsUsbNet ()
- !CMcsUsbNet ()
- DeviceEnumNet GetDeviceEnum ()
- virtual uint32_t Connect ()

Opens a connection to the device.

virtual uint32_t Connect (unsigned int LockMask)

Opens a connection to the device.

virtual uint32_t Connect (CMcsUsbListEntryNet[^] entry)

Opens a connection to the device.

virtual uint32_t Connect (CMcsUsbListEntryNet[^] entry, unsigned int LockMask)

Opens a connection to the device.

- virtual uint32 t GetStatus ([System::Runtime::InteropServices::Out]uint32 t% iStatus)
- virtual bool IsConnected ()

Check if a device is Connected.

virtual void Disconnect ()

Disconnect from a device.

- CMcsUsbListEntryNet ^ GetUsbListEntry ()
- virtual String \(^\) GetSerialNumber ()

Query the Serial Number of the device.

- DriverVersionNet [^] GetVersion ()
- DriverVersionNet ^ GetVersion (CFirmwareDestinationNet dest)
- DeviceIdNet ^ GetDeviceId ()
- uint32_t GetIdent ([System::Runtime::InteropServices::Out]String^% Answer)
- void MultibootSelectImage (unsigned int sector)

Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.

String \(^\) MultibootGetImageId (unsigned int sector)

Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.

uint32 t MultibootGetCypressImageId (unsigned int sector)

Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.

uint32_t MultibootGetSelectedImage ()

Gets sector index of selected FPGA boot image on IFB

uint32_t GetMea21UsbPort ()

Gets the USB port if an IFB that is used by this connection

HeadstageIdEnumNet GetHeadstageID (uint32_t headstage)

Gets the ID of a connected headstage.

bool GetHeadstagePresent (uint32 t headstage)

queries whether a headstage is present

bool GetHeadstageActive (uint32_t headstage)

queries whether a headstage is active

void RescanHeadstage (uint32 t headstage)

rescans and activates a headstage

- array< BYTE > ^ GetSoftwareKey (unsigned int index)
- void SetSoftwareKey (unsigned int index, array< BYTE >^ buffer)
- void RemoveSoftwareKey (unsigned int index)
- void AddSoftwareKey (String[^] key)
- bool EmptyKey (String[^] key)
- bool ValidKey (String[^] key, [System::Runtime::InteropServices::Out]String[^]% serial_number)
- bool HasSoftwareKey (uint8_t ProgrammID, uint8_t majorversion)
- bool HasSoftwareKey (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- String \(^{\text{GetSoftwareKeyString}}\) (uint8 t ProgrammID, uint8 t majorversion)
- String ^ GetSoftwareKeyString (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- bool IsDeviceHighSpeedCapable ()
- bool IsDeviceHighSpeed ()
- McsUsbSpeedEnumNet GetDeviceCapableSpeed ()
- McsUsbSpeedEnumNet GetDeviceSpeed ()

Query the Connection Speed of the device.

- unsigned int TxnTestMemoryWrite (unsigned short index)
- unsigned int TxnTestMemoryReadAndCheck (unsigned short index)
- void TxnSetSerialNumber (unsigned int number)
- unsigned int TxnGetSerialNumber ()
- unsigned int ReadRegister (unsigned int reg)

- array< uint32 t > ^ ReadRegister (unsigned int reg, int length)
- unsigned int ReadRegister32 (unsigned int adr)
- unsigned int ReadRegisterTimeSlot (unsigned int reg, int TimeSlot)
- void WriteRegister (unsigned int reg, unsigned int value)
- void WriteRegisterValue (unsigned int reg, unsigned int value)
- void WriteRegister32 (unsigned int adr, unsigned int value)
- void WriteRegister (unsigned int reg, array< unsigned int >[∧] values)
- void WriteRegisterArray (unsigned int reg, array< unsigned int >^ values)
- · void WriteRegisterTimeSlot (unsigned int reg, unsigned int value, int TimeSlot)
- void WriteRegisterTimeSlot (unsigned int reg, array< unsigned int >[^] values, int TimeSlot)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::← InteropServices::Out]uint32_t% DMA_value)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::
 —
 InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::
 —
 InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void WriteEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value)
- void WriteEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize)
- void WriteEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- unsigned int GetLastUSBError ()
- void ThrowCUsbExceptionNetOnError (uint32_t status)
- uint32_t GetDeviceCannotStallOutRequests ()
- String ^ GetHardwareRevision ()
- unsigned int GetFirmwareVersion (CFirmwareDestinationNet destination)

Gets the firmware version for the destination.

- uint8 t GetNumConfigurations ()
- uint8 t GetConfiguration ()
- void SetConfiguration (uint8 t config)
- uint32 t GetDeviceRootHubVendorID ()

Gets the Vendor ID of the USB root hub the device is connected to.

UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ()

Gets the Vendor ID of the USB root hub the device is connected to.

String \(^\) GetDeviceRootHubVendorName ()

Gets the Vendor Name of the USB root hub the device is connected to.

void EnableExceptions (bool enable)

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.

- bool IsExceptionsEnabled ()
- uint32_t GetStatusOfLastCommand ()

Gets the status of the last call to the McsUsb Library.

void AssociateToThis (CMcsUsbNet[^] device)

Static Public Member Functions

static String \(^\) GetErrorText (unsigned int Status)

Gets the error text string that belongs to a status number.

Static Public Attributes

```
    static const uint32 t Status Crc = (0xE0100001L)

    static const uint32 t Status Btstuff = (0xE0100002L)

    static const uint32 t Status DataToggleMismatch = (0xE0100003L)

    static const uint32 t Status Stall = (0xE0100004L)

    static const uint32_t Status_DevNotResponding = (0xE0100005L)

    static const uint32 t Status PidCheckFailure = (0xE0100006L)

    static const uint32_t Status_UnexpectedPid = (0xE0100007L)

• static const uint32_t Status_DataOverrun = (0xE0100008L)

    static const uint32 t Status DataUnderrun = (0xE0100009L)

    static const uint32 t Status BufferOverrun = (0xE010000CL)

    static const uint32 t Status BufferUnderrun = (0xE010000DL)

    static const uint32 t Status NotAccessed = (0xE010000FL)

    static const uint32_t Status_Fifo = (0xE0100010L)

    static const uint32 t Status EndpointHalted = (0xE0100030L)

    static const uint32 t Status NoMemory = (0xE0100100L)

    static const uint32_t Status_InvalidUrbFunction = (0xE0100200L)

    static const uint32 t Status InvalidParameter = (0xE0100300L)

    static const uint32_t Status_ErrorBusy = (0xE0100400L)

    static const uint32 t Status RequestFailed = (0xE0100500L)

    static const uint32 t Status InvalidPipeHandle = (0xE0100600L)

    static const uint32_t Status_NoBandwidth = (0xE0100700L)

    static const uint32 t Status InternalHcError = (0xE0100800L)

    static const uint32 t Status ErrorShortTransfer = (0xE0100900L)

    static const uint32_t Status_BadStartFrame = (0xE0100A00L)

    static const uint32 t Status IsochRequestFailed = (0xE0100B00L)

    static const uint32 t Status FrameControlOwned = (0xE0100C00L)

    static const uint32 t Status ControlNotOwned = (0xE0100D00L)

    static const uint32_t Status_Canceled = (0xE0110000L)

    static const uint32_t Status_Canceling = (0xE0120000L)

    static const uint32 t Status AlreadyConfigured = (0xE0110001L)

    static const uint32 t Status Unconfigured = (0xE0110002L)

    static const uint32 t Status NoSuchDevice = (0xE01F0002L)

    static const uint32 t Status DeviceNotFound = (0xE01F0003L)

    static const uint32_t Status_NotSupported = (0xE01F0005L)

    static const uint32 t Status loPending = (0xE01F0006L)

    static const uint32 t Status IoTimeout = (0xE01F0007L)

    static const uint32 t Status DeviceRemoved = (0xE01F0008L)

    static const uint32 t Status PipeNotLinked = (0xE01F0009L)

    static const uint32_t Status_ConnectedPipes = (0xE01F000AL)

    static const uint32_t Status_DeviceLocked = (0xE01F0010L)

    static const uint32 t Status RequestMutexTimeout = (0xE01F0020L)

    static const uint32_t Status_RequestMutexFailed = (0xE01F0021L)

    static const uint32 t Status LastUsbErrorMismatch = (0xE01F0022L)

    static const uint32 t WPAError ScanningIsPending = ( (0xA0220000L) | 0x0036 )
```

Properties

• virtual String Serial Number [get]

11.56.1 Detailed Description

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

11.56.2 Constructor & Destructor Documentation

```
11.56.2.1 CMcsUsbNet() [1/2] CMcsUsbNet ( )
```

Initializes a new instance of the base class to handle MCS USB devices.

```
11.56.2.2 CMcsUsbNet() [2/2] CMcsUsbNet (

McsBusTypeEnumNet bustype)
```

Initializes a new instance of the base class to handle MCS USB devices.

Parameters

```
bustype Type of device to use, either USB or PCI.
```

```
11.56.2.3 \sim CMcsUsbNet() virtual \sim CMcsUsbNet () [virtual]
```

```
11.56.2.4 "!CMcsUsbNet() !CMcsUsbNet ( )
```

11.56.3 Member Function Documentation

```
11.56.3.1 AddSoftwareKey() void AddSoftwareKey ( String^{\land} key)
```

```
11.56.3.2 AssociateToThis() void AssociateToThis (

CMcsUsbNet^ device )
```

11.56.3.3 Connect() [1/4] virtual uint32_t Connect () [virtual]

Opens a connection to the device.

Returns

Error Status. 0 on success.

```
11.56.3.4 Connect() [2/4] virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry ) [virtual]
```

Opens a connection to the device.

Parameters

	entry	The Device List Entry for the device to be connected.
--	-------	---

Returns

Error Status. 0 on success.

```
11.56.3.5 Connect() [3/4] virtual uint32_t Connect (

CMcsUsbListEntryNet^ entry,

unsigned int LockMask) [virtual]
```

Opens a connection to the device.

Parameters

entry	The Device List Entry for the device to be connected.
LockMask	The Lock Mask for this connection.

Returns

Error Status. 0 on success.

```
11.56.3.6 Connect() [4/4] virtual uint32_t Connect (
unsigned int LockMask) [virtual]
```

Opens a connection to the device.

Parameters

LockMask The Lock Mask for this connection	۱.
--	----

Returns

Error Status. 0 on success.

```
11.56.3.7 Disconnect() virtual void Disconnect ( ) [virtual]
```

Disconnect from a device.

```
11.56.3.8 EmptyKey() bool EmptyKey ( String^{\wedge} key )
```

```
11.56.3.9 EnableExceptions() void EnableExceptions ( bool enable)
```

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.

Parameters

```
enable True to enable Exceptions, False to disable.
```

```
11.56.3.10 EraseEepromRegisterPreconfig() [1/3] void EraseEepromRegisterPreconfig ( uint32_t EEPROMBase,
```

```
uint32_t DMA_reg )
```

```
11.56.3.11 EraseEepromRegisterPreconfig() [2/3] void EraseEepromRegisterPreconfig (
```

```
uint32_t EEPROMBase,
uint32_t DMA_reg,
uint32_t EEPROMSize )
```

```
\textbf{11.56.3.12} \quad \textbf{EraseEepromRegisterPreconfig() [3/3]} \quad \texttt{void EraseEepromRegisterPreconfig ()} \\
               uint32_t EEPROMBase,
               uint32_t DMA_reg,
               uint32_t EEPROMSize,
               uint32_t EepromStartAddress )
11.56.3.13 GetConfiguration() uint8_t GetConfiguration ( )
\textbf{11.56.3.14} \quad \textbf{GetDeviceCannotStallOutRequests()} \quad \texttt{uint32\_t} \quad \texttt{GetDeviceCannotStallOutRequests} \quad \textbf{()}
11.56.3.15 GetDeviceCapableSpeed() McsUsbSpeedEnumNet GetDeviceCapableSpeed ( )
11.56.3.16 GetDeviceEnum() DeviceEnumNet GetDeviceEnum ( )
11.56.3.17 GetDeviceId() DeviceIdNet ^ GetDeviceId ( )
11.56.3.18 GetDeviceRootHubVendorEnum() UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ( )
Gets the Vendor ID of the USB root hub the device is connected to.
Returns
     An enum which enumerates the PCI Vendor ID.
11.56.3.19 GetDeviceRootHubVendorID() uint32_t GetDeviceRootHubVendorID ( )
Gets the Vendor ID of the USB root hub the device is connected to.
Returns
```

The PCI Vendor ID, 0x8086 for Intel, 0x1912 for Renesas, 0x1b21 for ASMedia.

11.56.3.20 GetDeviceRootHubVendorName() String ^ GetDeviceRootHubVendorName ()

Gets the Vendor Name of the USB root hub the device is connected to.

Returns

The PCI Vendor Name, either "Intel", "Renesas", "ASMedia" or "unknown".

11.56.3.21 GetDeviceSpeed() McsUsbSpeedEnumNet GetDeviceSpeed ()

Query the Connection Speed of the device.

Returns

0 for Low-Speed, 1 for Full-Speed, 2 for High-Speed and 3 for SuperSpeed.

Gets the error text string that belongs to a status number.

Parameters

Status	The status number you want the text for.
--------	--

Returns

The error text string that belongs to the status number.

11.56.3.23 GetFirmwareVersion() unsigned int GetFirmwareVersion (CFirmwareDestinationNet *destination*)

Gets the firmware version for the destination.

Parameters

destination	η.	The destination to be queried.

Returns

The firmware version as a 32 bit number, the upper 16 bit contain the majaor version number, the lower 16 bit the minor version number.

11.56.3.24 GetHardwareRevision() String $^{\wedge}$ GetHardwareRevision ()

```
11.56.3.25 GetHeadstageActive() bool GetHeadstageActive ( uint32_t headstage )
```

queries whether a headstage is active

Parameters

in headstage the headstage nu	umber (0 or 1)
-------------------------------	----------------

Returns

true if the headstage is active

11.56.3.26 GetHeadstagelD() HeadstageIdEnumNet GetHeadstageID (uint32_t headstage)

Gets the ID of a connected headstage.

Parameters

ſ			
١	in	headstage	the headstage number (0 or 1)

Returns

enumerated Headstage ID

11.56.3.27 GetHeadstagePresent() bool GetHeadstagePresent (uint32_t headstage)

queries whether a headstage is present

Parameters

in	headstage	the headstage number (0 or 1)

Returns

true if the headstage is present

```
11.56.3.28 GetIdent() uint32_t GetIdent (
              [System::Runtime::InteropServices::Out] String^% Answer )
11.56.3.29 GetLastUSBError() unsigned int GetLastUSBError ()
11.56.3.30 GetMea21UsbPort() uint32_t GetMea21UsbPort ( )
Gets the USB port if an IFB that is used by this connection
Returns
     number of used port; range: 0..1
11.56.3.31 GetNumConfigurations() uint8_t GetNumConfigurations ( )
11.56.3.32 GetSerialNumber() virtual String ^ GetSerialNumber ( ) [virtual]
Query the Serial Number of the device.
Returns
     The Serial Number.
11.56.3.33 GetSoftwareKey() array<br/>SYTE> ^{\wedge} GetSoftwareKey (
             unsigned int index)
11.56.3.34 GetSoftwareKeyString() [1/2] String ^ GetSoftwareKeyString (
             SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
             uint8_t majorversion )
11.56.3.35 GetSoftwareKeyString() [2/2] String ^ GetSoftwareKeyString (
             uint8_t ProgrammID,
             uint8_t majorversion )
```

```
11.56.3.36 GetStatus() virtual uint32_t GetStatus (
             [System::Runtime::InteropServices::Out] uint32_t% iStatus ) [virtual]
11.56.3.37 GetStatusOfLastCommand() uint32_t GetStatusOfLastCommand ( )
Gets the status of the last call to the McsUsb Library.
Returns
    The Error Status of the last McsUsb command. 0 on success.
11.56.3.38 GetUsbListEntry() CMcsUsbListEntryNet ^ GetUsbListEntry ( )
11.56.3.39 GetVersion() [1/2] DriverVersionNet ^ GetVersion ( )
11.56.3.40 GetVersion() [2/2] DriverVersionNet ^ GetVersion (
             CFirmwareDestinationNet dest )
11.56.3.41 HasSoftwareKey() [1/2] bool HasSoftwareKey (
             SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
             uint8_t majorversion )
11.56.3.42 HasSoftwareKey() [2/2] bool HasSoftwareKey (
             uint8_t ProgrammID,
             uint8_t majorversion )
11.56.3.43 IsConnected() virtual bool IsConnected ( ) [virtual]
Check if a device is Connected.
```

Generated by Doxygen

true if the device is connected.

Returns

```
11.56.3.44 IsDeviceHighSpeed() bool IsDeviceHighSpeed ( )
```

```
11.56.3.45 IsDeviceHighSpeedCapable() bool IsDeviceHighSpeedCapable ( )
```

```
11.56.3.47 MultibootGetCypressImageId() uint32_t MultibootGetCypressImageId ( unsigned int sector )
```

Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.

Returns

The magic ident code of the image.

```
11.56.3.48 MultibootGetImageId() String ^ MultibootGetImageId ( unsigned int sector )
```

Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.

Returns

The magic ident code of the image.

```
11.56.3.49 MultibootGetSelectedImage() uint32_t MultibootGetSelectedImage ( )
```

Gets sector index of selected FPGA boot image on IFB

Returns

Sector index of image; range: 0..2

```
11.56.3.50 MultibootSelectImage() void MultibootSelectImage ( unsigned int sector )
```

Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.

Returns

```
Throws exception on error.
\textbf{11.56.3.51} \quad \textbf{ReadEepromRegisterPreconfig() [1/3]} \quad \texttt{bool ReadEepromRegisterPreconfig ()} \quad \textbf{[1/3]} \quad \textbf{
                                                            uint32_t EEPROMBase,
                                                            uint32_t DMA_reg,
                                                             [System::Runtime::InteropServices::Out] uint32_t% DMA_value )
11.56.3.52 ReadEepromRegisterPreconfig() [2/3] bool ReadEepromRegisterPreconfig (
                                                           uint32_t EEPROMBase,
                                                            uint32_t DMA_reg,
                                                             [System::Runtime::InteropServices::Out] uint32_t% DMA_value,
                                                            uint32_t EEPROMSize )
11.56.3.53 ReadEepromRegisterPreconfig() [3/3] bool ReadEepromRegisterPreconfig (
                                                            uint32_t EEPROMBase,
                                                            uint32_t DMA_reg,
                                                            [System::Runtime::InteropServices::Out] uint32_t% DMA_value,
                                                            uint32_t EEPROMSize,
                                                            uint32_t EepromStartAddress )
11.56.3.54 ReadRegister() [1/2] unsigned int ReadRegister (
                                                            unsigned int reg )
11.56.3.55 ReadRegister() [2/2] array<uint32_t> ^ ReadRegister (
                                                            unsigned int reg,
                                                            int length )
11.56.3.56 ReadRegister32() unsigned int ReadRegister32 (
```

unsigned int adr)

```
11.56.3.57 ReadRegisterTimeSlot() unsigned int ReadRegisterTimeSlot (
              unsigned int reg,
              int TimeSlot )
11.56.3.58 RemoveSoftwareKey() void RemoveSoftwareKey (
             unsigned int index )
\textbf{11.56.3.59} \quad \textbf{RescanHeadstage()} \quad \texttt{void RescanHeadstage ()}
             uint32_t headstage )
rescans and activates a headstage
Parameters
       headstage
                    the headstage number (0 or 1)
   in
11.56.3.60 SetConfiguration() void SetConfiguration (
             uint8_t config )
11.56.3.61 SetSoftwareKey() void SetSoftwareKey (
             unsigned int index,
              array < BYTE >^{\wedge} buffer )
11.56.3.62 ThrowCUsbExceptionNetOnError() void ThrowCUsbExceptionNetOnError (
             uint32_t status )
11.56.3.63 TxnGetSerialNumber() unsigned int TxnGetSerialNumber ()
11.56.3.64 TxnSetSerialNumber() void TxnSetSerialNumber (
              unsigned int number )
```

```
11.56.3.65 TxnTestMemoryReadAndCheck() unsigned int TxnTestMemoryReadAndCheck (
             unsigned short index)
11.56.3.66 TxnTestMemoryWrite() unsigned int TxnTestMemoryWrite (
             unsigned short index)
11.56.3.67 ValidKey() [1/2] bool ValidKey (
             String^{\wedge} key,
             [System::Runtime::InteropServices::Out] \ String^{\ } \ serial\_number \ )
11.56.3.68 ValidKey() [2/2] bool ValidKey (
             String^{\wedge} key,
             uint8_t ProgrammID,
             uint8_t majorversion,
             [System::Runtime::InteropServices::Out] String^% serial_number )
11.56.3.69 WriteEepromRegisterPreconfig() [1/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value )
11.56.3.70 WriteEepromRegisterPreconfig() [2/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value,
             uint32_t EEPROMSize )
11.56.3.71 WriteEepromRegisterPreconfig() [3/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value,
             uint32_t EEPROMSize,
             uint32_t EepromStartAddress )
```

```
11.56.3.72 WriteRegister() [1/2] void WriteRegister (
             unsigned int reg,
             array< unsigned int >^{\wedge} values)
11.56.3.73 WriteRegister() [2/2] void WriteRegister (
             unsigned int reg,
             unsigned int value )
11.56.3.74 WriteRegister32() void WriteRegister32 (
             unsigned int adr,
             unsigned int value )
11.56.3.75 WriteRegisterArray() void WriteRegisterArray (
             unsigned int reg,
             array< unsigned int >^{\land} values)
11.56.3.76 WriteRegisterTimeSlot() [1/2] void WriteRegisterTimeSlot (
             unsigned int reg,
             array< unsigned int >^{\wedge} values,
             int TimeSlot )
11.56.3.77 WriteRegisterTimeSlot() [2/2] void WriteRegisterTimeSlot (
             unsigned int reg,
             unsigned int value,
             int TimeSlot )
11.56.3.78 WriteRegisterValue() void WriteRegisterValue (
             unsigned int reg,
             unsigned int value )
```

11.56.4 Member Data Documentation

```
11.56.4.1 Status_AlreadyConfigured const uint32_t Status_AlreadyConfigured = (0xE0110001L)
[static]
11.56.4.2 Status_BadStartFrame const uint32_t Status_BadStartFrame = (0xE0100A00L) [static]
11.56.4.3 Status_Btstuff const uint32_t Status_Btstuff = (0xE0100002L) [static]
11.56.4.4 Status_BufferOverrun const uint32_t Status_BufferOverrun = (0xE010000CL) [static]
11.56.4.5 Status_BufferUnderrun const uint32_t Status_BufferUnderrun = (0xE010000DL) [static]
11.56.4.6 Status_Canceled const uint32_t Status_Canceled = (0xE0110000L) [static]
11.56.4.7 Status_Canceling const uint32_t Status_Canceling = (0xE0120000L) [static]
11.56.4.8 Status_ConnectedPipes const uint32_t Status_ConnectedPipes = (0xE01F000AL) [static]
```

11.56.4.11 Status_DataOverrun const uint32_t Status_DataOverrun = (0xE0100008L) [static]

11.56.4.10 Status_Crc const uint32_t Status_Crc = (0xE01000001L) [static]

11.56.4.9 Status_ControlNotOwned const uint32_t Status_ControlNotOwned = (0xE0100D00L) [static]

```
11.56.4.12 Status_DataToggleMismatch const uint32_t Status_DataToggleMismatch = (0xE0100003L)
[static]
11.56.4.13 Status_DataUnderrun const uint32_t Status_DataUnderrun = (0xE0100009L) [static]
11.56.4.14 Status_DeviceLocked const uint32_t Status_DeviceLocked = (0xE01F0010L) [static]
11.56.4.15 Status_DeviceNotFound const uint32_t Status_DeviceNotFound = (0xE01F0003L) [static]
11.56.4.16 Status DeviceRemoved const uint32_t Status_DeviceRemoved = (0xE01F0008L) [static]
11.56.4.17 Status_DevNotResponding const uint32_t Status_DevNotResponding = (0xE0100005L)
[static]
11.56.4.18 Status_EndpointHalted const uint32_t Status_EndpointHalted = (0xE0100030L) [static]
11.56.4.19 Status_ErrorBusy const uint32_t Status_ErrorBusy = (0xE0100400L) [static]
11.56.4.20 Status_ErrorShortTransfer const uint32_t Status_ErrorShortTransfer = (0xE0100900L)
[static]
11.56.4.21 Status_Fifo const uint32_t Status_Fifo = (0xE0100010L) [static]
```

11.56.4.22 Status_FrameControlOwned const uint32_t Status_FrameControlOwned = (0xE0100C00L)

[static]

```
11.56.4.23 Status_InternalHcError const uint32_t Status_InternalHcError = (0xE0100800L) [static]
11.56.4.24 Status_InvalidParameter const uint32_t Status_InvalidParameter = (0xE0100300L) [static]
11.56.4.25 Status_InvalidPipeHandle const uint32_t Status_InvalidPipeHandle = (0xE0100600L)
 [static]
11.56.4.26 Status_InvalidUrbFunction const uint32_t Status_InvalidUrbFunction = (0xE0100200L)
 [static]
11.56.4.27 Status_loPending const uint32_t Status_IoPending = (0xE01F0006L) [static]
11.56.4.28 Status_IoTimeout const uint32_t Status_IoTimeout = (0xE01F0007L) [static]
11.56.4.29 Status_IsochRequestFailed const uint32_t Status_IsochRequestFailed = (0xE0100B00L)
 [static]
\textbf{11.56.4.30} \quad \textbf{Status\_LastUsbErrorMismatch} \quad \texttt{const uint} \\ 32\_t \quad \texttt{Status\_LastUsbErrorMismatch} \\ = (0xE01 \leftarrow 0xE01) \\ + (0xE01 \leftarrow 0xE01) \\
F0022L) [static]
11.56.4.31 Status NoBandwidth const uint32_t Status_NoBandwidth = (0xE0100700L) [static]
11.56.4.32 Status_NoMemory const uint32_t Status_NoMemory = (0xE0100100L) [static]
11.56.4.33 Status_NoSuchDevice const uint32_t Status_NoSuchDevice = (0xE01F0002L) [static]
```

```
11.56.4.34 Status_NotAccessed const uint32_t Status_NotAccessed = (0xE010000FL) [static]
11.56.4.35 Status_NotSupported const uint32_t Status_NotSupported = (0xE01F0005L) [static]
11.56.4.36 Status_PidCheckFailure const uint32_t Status_PidCheckFailure = (0xE0100006L) [static]
11.56.4.37 Status_PipeNotLinked const uint32_t Status_PipeNotLinked = (0xE01F0009L) [static]
11.56.4.38 Status RequestFailed const uint32_t Status_RequestFailed = (0xE0100500L) [static]
11.56.4.39 Status RequestMutexFailed const uint32_t Status_RequestMutexFailed = (0xE01F0021L)
[static]
11.56.4.40 Status_RequestMutexTimeout const uint32_t Status_RequestMutexTimeout = (0xE01↔
F0020L) [static]
11.56.4.41 Status_Stall const uint32_t Status_Stall = (0xE0100004L) [static]
11.56.4.42 Status Unconfigured const uint32_t Status_Unconfigured = (0xE0110002L) [static]
11.56.4.43 Status_UnexpectedPid const uint32_t Status_UnexpectedPid = (0xE0100007L) [static]
11.56.4.44 WPAError ScanningIsPending const uint32_t WPAError_ScanningIsPending = ( (0x↔
```

A0220000L) | 0x0036) [static]

11.56.5 Property Documentation


11.56.5.1 SerialNumber virtual String SerialNumber [get]

11.57 CMcsUsbPointerContainer Class Reference

11.58 CMEA2100x256FunctionNet Class Reference

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

Inheritance diagram for CMEA2100x256FunctionNet:

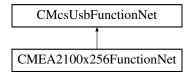

Public Member Functions

CMEA2100x256FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pME
 — A2100x256FunctionPointerContainer)

Initializes a new instance of the CMEA2100x256FunctionNet class.

- CMEA2100x256FunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CMEA2100x256FunctionNet ()
- !CMEA2100x256FunctionNet ()
- StimulationLayoutConfigurationEnumNet GetLayoutConfiguration ()

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of $D \leftarrow AC$ channels available per well is Mcs::Usb::Cstg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::Cstg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

• void SetLayoutConfiguration (StimulationLayoutConfigurationEnumNet LayoutConfiguration)

Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D← AC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

Additional Inherited Members

11.58.1 Detailed Description

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

11.58.2 Constructor & Destructor Documentation

```
11.58.2.1 CMEA2100x256FunctionNet() [1/2] CMEA2100x256FunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pMEA2100x256FunctionPointerContainer)
```

Initializes a new instance of the CMEA2100x256FunctionNet class.

```
11.58.2.2 CMEA2100x256FunctionNet() [2/2] CMEA2100x256FunctionNet (
CMcsUsbNet^ mcsusb )
```

11.58.2.3 ~CMEA2100x256FunctionNet() virtual ~CMEA2100x256FunctionNet () [virtual]

```
11.58.2.4 "!CMEA2100x256FunctionNet() !CMEA2100x256FunctionNet ()
```

11.58.3 Member Function Documentation

```
11.58.3.1 GetLayoutConfiguration() StimulationLayoutConfigurationEnumNet GetLayoutConfiguration ()
```

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D \leftarrow AC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

Returns

The currently active stimulation layout configuration.

```
11.58.3.2 SetLayoutConfiguration() void SetLayoutConfiguration (
StimulationLayoutConfigurationEnumNet LayoutConfiguration)
```


Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of D \leftarrow AC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

Parameters

Laurant Cambian matian	The new attraction levels acostic matter
LavouiConilguration	The new stimulation layout configuration.

11.59 CMeaAudioFunctionNet Class Reference

Inheritance diagram for CMeaAudioFunctionNet:

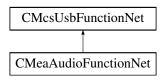

Classes

· struct s_setaudionet

Public Member Functions

- CMeaAudioFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] meaAudioFunction←
 PointerContainer)
- CMeaAudioFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual uint32 t GetNumberOfAudioChannels ()

Gets the number of available audio channels.

virtual uint32_t SetAudioChannels (array< s_setaudionet^>^ channels)

Sets the electrode to monitor and amplification for the audio channels.

- virtual uint32_t SetAudioChannels (array < s_setaudionet^>^ channels, unsigned int virtualDevice)
 - Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out]array< s_setaudionet^>^% channels)

Gets the electrode to monitor and amplification for the audio channels.

virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out]array< s_setaudionet^>^% channels, unsigned int virtualDevice)

Gets the electrode to monitor and amplification for the audio channels.

Additional Inherited Members

11.59.1 Constructor & Destructor Documentation

```
11.59.1.1 CMeaAudioFunctionNet() [1/2] CMeaAudioFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ meaAudioFunctionPointerContainer)
```

```
11.59.1.2 CMeaAudioFunctionNet() [2/2] CMeaAudioFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.59.2 Member Function Documentation

Gets the electrode to monitor and amplification for the audio channels.

Parameters

channels	Struct which contains the electrode (channel) and amplification on return.
----------	--

Returns

Error Status. 0 on success.

Gets the electrode to monitor and amplification for the audio channels.

Parameters

channels Struct which contains the electrode (channel) and ampli	lification on return.
--	-----------------------

Parameters

virtualDevice	Virtual device to use.
---------------	------------------------

Returns

Error Status. 0 on success.

11.59.2.3 GetNumberOfAudioChannels() virtual uint32_t GetNumberOfAudioChannels () [virtual]

Gets the number of available audio channels.

Returns

The number of audio channels available, 0 when there are none.

Sets the electrode to monitor and amplification for the audio channels.

Parameters

Returns

Error Status. 0 on success.

Sets the electrode to monitor and amplification for the audio channels.

Parameters

channels	Struct which defines the electrode (channel) and amplification.
----------	---

Parameters


Returns

Error Status. 0 on success.

11.60 CMeaCleanDeviceNet Class Reference

CMeaCleanDeviceNet is the class to access the MEA Clean device.

Inheritance diagram for CMeaCleanDeviceNet:

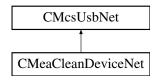

Public Member Functions

CMeaCleanDeviceNet ()

Initializes a new instance of the CMeaCleanDeviceNet class.

- virtual ~CMeaCleanDeviceNet ()
- !CMeaCleanDeviceNet ()
- void Start ()

Starts a MEA Clean run.

• void Stop ()

Stops a MEA Clean run.

void SetSlope (uint32_t voltageSlope)

Sets the voltage slope.

void SetCycles (uint32_t cycles)

Sets the number of cycles.

void SetMinVoltage (int32_t voltageMin)

Sets the lower voltage level.

void SetMaxVoltage (int32 t voltageMax)

Sets the upper voltage level.

• bool IsRunning ()

Gets if the MEA Clean device is running.

• uint32_t GetSlope ()

Gets the voltage slope.

uint32_t GetCycles ()

Gets the number of cycles.

• int32_t GetMinVoltage ()

Gets the lower voltage level.

int32_t GetMaxVoltage ()

Gets the upper voltage level

• int32_t GetOutputVoltage ()

Gets the output voltage.

int32_t GetCycle ()
 Gets the current cycle.

Additional Inherited Members

11.60.1 Detailed Description

CMeaCleanDeviceNet is the class to access the MEA Clean device.

11.60.2 Constructor & Destructor Documentation

11.60.2.1 CMeaCleanDeviceNet() CMeaCleanDeviceNet ()

Initializes a new instance of the CMeaCleanDeviceNet class.

```
11.60.2.2 ~CMeaCleanDeviceNet() virtual ~CMeaCleanDeviceNet () [virtual]
11.60.2.3 "!CMeaCleanDeviceNet() !CMeaCleanDeviceNet ()
11.60.3 Member Function Documentation
11.60.3.1 GetCycle() int32_t GetCycle ()
Gets the current cycle.
Returns
     The cycle number.
11.60.3.2 GetCycles() uint32_t GetCycles ()
Gets the number of cycles.
Returns
     The number of cycles to run for.
11.60.3.3 GetMaxVoltage() int32_t GetMaxVoltage ( )
Gets the upper voltage level
Returns
     The upper voltage level in mV.
11.60.3.4 GetMinVoltage() int32_t GetMinVoltage ()
Gets the lower voltage level.
Returns
```

The lower voltage level in mV.

11.60.3.5 GetOutputVoltage() int32_t GetOutputVoltage ()

Gets the output voltage.

Returns

The output voltage in mV.

```
11.60.3.6 GetSlope() uint32_t GetSlope ()
```

Gets the voltage slope.

Returns

The voltage slope in mV/s.

11.60.3.7 IsRunning() bool IsRunning ()

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

```
11.60.3.8 SetCycles() void SetCycles (
uint32_t cycles)
```

Sets the number of cycles.

Parameters

cycles The number of cycles to run for (0 .. 99).

Sets the upper voltage level.

Parameters

voltaneMay	The upper voltage level in mV ((-16 16 V/)
vuitautiviax	THE UDDEL VUILAGE IEVEL III III V (-1.0 1.0 v /.

Sets the lower voltage level.

Parameters

voltageMin The lower voltage level in mV (-1.6 .. 1.6 V).

```
11.60.3.11 SetSlope() void SetSlope (
uint32_t voltageSlope )
```

Sets the voltage slope.

Parameters

voltageSlope	The voltage slope in mV/s (range 0 60 V/s).
--------------	---

```
11.60.3.12 Start() void Start ()
```

Starts a MEA Clean run.


```
11.60.3.13 Stop() void Stop ()
```

Stops a MEA Clean run.

11.61 CMeaCoatDeviceNet Class Reference

CMeaCoatDeviceNet is the class to access the MEA Coat device.

Inheritance diagram for CMeaCoatDeviceNet:

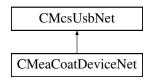

Public Member Functions

CMeaCoatDeviceNet ()

Initializes a new instance of the CMeaCoatDeviceNet class.

- virtual ∼CMeaCoatDeviceNet ()
- !CMeaCoatDeviceNet ()
- void Start ()

Starts a MEA Coat run.

• void Stop ()

Stops a MEA Coat run.

void SetSlope (int32_t currentSlope)

Sets the current slope.

void SetDuration (uint32_t duration)

Sets the duration of a MEA Coat run.

void SetMaxCurrent (uint32_t currentMax)

Sets the limit of the current ramp (absolute value).

void SetOffsetCurrent (int32 t currentOffset)

Sets the offset of the current.

• bool IsRunning ()

Gets if the MEA Clean device is running.

int32 t GetSlope ()

Gets the current slope.

uint32_t GetDuration ()

Gets the duration of a MEA Coat run.

uint32_t GetMaxCurrent ()

Gets the limit of the current ramp (absolute value).

• int32_t GetOffsetCurrent ()

Gets the offset of the current.

int32_t GetOutputCurrent ()

Gets the output current.

• int32_t GetTimeInPlateau ()

Gets the time in the plateau.

• void SetPauseDuration (uint32 t pauseDuration)

Sets the duration of the pause between MEA Coat pulses.

uint32_t GetPauseDuration ()

Gets the duration of the pause between MEA Coat pulses.

• int32_t GetTimeInPause ()

Gets the time in the pause.

void SetCycles (uint32_t cycles)

Sets the number of cycles.

• uint32_t GetCycles ()

Gets the number of cycles.

int32_t GetCurrentCycle ()

Gets the current cycle.

Additional Inherited Members

11.61.1 Detailed Description

CMeaCoatDeviceNet is the class to access the MEA Coat device.

11.61.2 Constructor & Destructor Documentation

11.61.2.1 CMeaCoatDeviceNet() CMeaCoatDeviceNet () Initializes a new instance of the CMeaCoatDeviceNet class. 11.61.2.2 ~CMeaCoatDeviceNet() virtual ~CMeaCoatDeviceNet () [virtual] 11.61.2.3 "!CMeaCoatDeviceNet() !CMeaCoatDeviceNet () 11.61.3 Member Function Documentation 11.61.3.1 GetCurrentCycle() int32_t GetCurrentCycle () Gets the current cycle.

dets the current t

Returns

The cycle number.

```
11.61.3.2 GetCycles() uint32_t GetCycles ( )
```

Gets the number of cycles.

Returns

The number of cycles to run for.

11.61.3.3 GetDuration() uint32_t GetDuration ()

Gets the duration of a MEA Coat run.

Returns

The duration in ms.

11.61.3.4 GetMaxCurrent() uint32_t GetMaxCurrent ()

Gets the limit of the current ramp (absolute value).

Returns

The limit of the current ramp in pA (absolute value).

11.61.3.5 GetOffsetCurrent() int32_t GetOffsetCurrent ()

Gets the offset of the current.

Returns

The offset of the current in pA.

11.61.3.6 GetOutputCurrent() int32_t GetOutputCurrent ()

Gets the output current.

Returns

The output current in pA.

11.61.3.7 GetPauseDuration() uint32_t GetPauseDuration ()

Gets the duration of the pause between MEA Coat pulses.

Returns

The duration in ms.

11.61.3.8 GetSlope() int32_t GetSlope ()

Gets the current slope.

Returns

The current slope in pA/s.

11.61.3.9 GetTimeInPause() int32_t GetTimeInPause ()

Gets the time in the pause.

Returns

The time in the pause in ms.

11.61.3.10 GetTimeInPlateau() int32_t GetTimeInPlateau ()

Gets the time in the plateau.

Returns

The time in the plateau in ms.

11.61.3.11 IsRunning() bool IsRunning ()

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

```
11.61.3.12 SetCycles() void SetCycles (
uint32_t cycles)
```

Sets the number of cycles.

Parameters

cycles The number of cycles to run for (0 .. 99).

```
11.61.3.13 SetDuration() void SetDuration ( uint32_t duration)
```

Sets the duration of a MEA Coat run.

Parameters

on in ms (range 0 65 s).

```
11.61.3.14 SetMaxCurrent() void SetMaxCurrent ( uint32_t currentMax )
```

Sets the limit of the current ramp (absolute value).

Parameters

currentMax | The limit of the current ramp in pA (absolute value, 0 .. 18 nA).

Sets the offset of the current.

Parameters

```
currentOffset The offset of the current in pA (-10 .. 10 nA).
```

```
11.61.3.16 SetPauseDuration() void SetPauseDuration ( uint32_t pauseDuration )
```

Sets the duration of the pause between MEA Coat pulses.

Parameters

pauseDuration The duration in ms (range 0 .. 65 s).

```
11.61.3.17 SetSlope() void SetSlope ( int32_t currentSlope )
```

Sets the current slope.

Parameters

currentSlope The current slope in pA/s (range -65 .. 65 nA/s).

11.61.3.18 Start() void Start ()

Starts a MEA Coat run.


11.61.3.19 Stop() void Stop ()

Stops a MEA Coat run.

11.62 CMeaDeviceNet Class Reference

Base class for MEA data acquisition devices.

Inheritance diagram for CMeaDeviceNet:

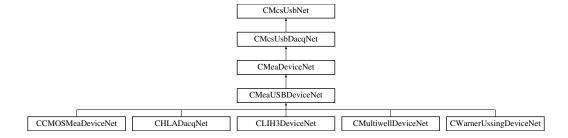

Public Member Functions

CMeaDeviceNet (McsBusTypeEnumNet bustype)

Initializes a new instance of CMeaDeviceNet class.

- CMeaDeviceNet (McsBusTypeEnumNet bustype, OnChannelData[^] channelData, OnError[^] error)
 - Initializes a new instance of CMeaDeviceNet class.
- ∼CMeaDeviceNet ()
- virtual int32_t GetGain ()

Gets the amplifier gain of the device.

- int32_t GetEnumerationSpeed ()
- virtual int32 t GetAnalogGain ()

Gets the gain of the analog inputs of the device.

• virtual uint32_t EnableDigitalIn (bool enable, unsigned int virtualDevice)

Enable the digital data word in the datastream.

- virtual uint32_t EnableDigitalIn (DigitalDatastreamEnableEnumNet enable, unsigned int virtualDevice)
 - Enable digital data words in the datastream.
- virtual uint32_t EnableTimestamp (bool enable, unsigned int virtualDevice)

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

virtual uint32 t EnableChecksum (bool enable, unsigned int virtualDevice)

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

virtual void SetDigitalOut (unsigned int digout_value, int pulselength)

Generate a pulse on the digital output.

• virtual uint32_t SetNumberOfChannels (int NumberOfChannels)

Sets the number of analog channels in the datastream.

• virtual uint32_t SetNumberOfChannels (int NumberOfChannels, unsigned int virtualDevice)

Sets the number of analog channels in the datastream.

virtual uint32_t SetNumberOfAnalogChannels (unsigned int NumberOfChannels_HS1, unsigned int NumberOfChannels_HS2, unsigned int NumberOfChannels_DSP, unsigned int NumberOfChannels_IF, unsigned int virtualDevice)

Sets the number of analog channels in the datastream for the MEA2100 device.

- virtual uint32_t SetTriggerPeriod (int samples, unsigned int virtualDevice)
 - Sets the maximum number of samples per trigger.
- virtual uint32_t SetTriggerMaskValue (unsigned int mask, unsigned int value, unsigned int virtualDevice)

Defines a pattern on the digital dataword which will start a trigger when found.

Properties

- CWClassicFunctionNet[^] WClassicFunctionNet [get]
- CW2100 FunctionNet[^] W2100 FunctionNet [get]
- CMeaAudioFunctionNet[^] MeaAudioFunctionNet [get]
- CMeaDigitalDataFunctionNet[^] MeaDigitalDataFunctionNet [get]
- CMeaFeedbackFunctionNet[^] MeaFeedbackFunctionNet [get]
- virtual int Gain [get]

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

virtual int AnalogGain [get]

The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Additional Inherited Members

11.62.1 Detailed Description

Base class for MEA data acquisition devices.

There are two different device types for MEA data aquistion devices. There are the USB-MEA devices and the MC← _Card. In .NET both classes can be accessed by the contructor of the base class CMeaDeviceNet, which contructs the correct underlying C++ class for the USB-MEA device on the one hand or the MC_Card device on the other hand. Through this interface both device types USB-MEA devices and MC_Card devices can be accessed

11.62.2 Constructor & Destructor Documentation

```
11.62.2.1 CMeaDeviceNet() [1/2] CMeaDeviceNet (
McsBusTypeEnumNet bustype)
```

Initializes a new instance of CMeaDeviceNet class.

use, either USB or PCI.	bustype
-------------------------	---------

11.62.2.2 CMeaDeviceNet() [2/2] CMeaDeviceNet (

```
McsBusTypeEnumNet bustype,
OnChannelData^ channelData,
OnError^ error )
```

Initializes a new instance of CMeaDeviceNet class.

Parameters

bustype	Type of device to use, either USB or PCI.
---------	---

Parameters

Parameters

```
error | Callback to call when an error occurred.
```

11.62.2.3 \sim CMeaDeviceNet() \sim CMeaDeviceNet ()

11.62.3 Member Function Documentation

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable the digital data word in the datastream.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable digital data words in the datastream.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

11.62.3.5 GetAnalogGain() virtual int32_t GetAnalogGain () [virtual]

Gets the gain of the analog inputs of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.62.3.6 GetEnumerationSpeed() int32_t GetEnumerationSpeed ()

```
11.62.3.7 GetGain() virtual int32_t GetGain ( ) [virtual]
```

Gets the amplifier gain of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Generate a pulse on the digital output.

Parameters

pulselength	Pulselength in ms.
-------------	--------------------

11.62.3.9 SetNumberOfAnalogChannels() virtual uint32_t SetNumberOfAnalogChannels (

```
unsigned int NumberOfChannels_HS1,
unsigned int NumberOfChannels_HS2,
unsigned int NumberOfChannels_DSP,
unsigned int NumberOfChannels_IF,
unsigned int virtualDevice ) [virtual]
```

Sets the number of analog channels in the datastream for the MEA2100 device.

Parameters

NumberOfChannels_HS1	Number of analog channels from the Headstage 1.
----------------------	---

Parameters

NumberOfChannels_HS2	Number of analog channels from the Headstage 2.
----------------------	---

Parameters

NumberOfChannels_DSP	Number of data words from the DSP.
----------------------	------------------------------------

Parameters

NumberOfChannels⇔	Number of analog channels from the Interfaceboard.
_IF	

Parameters

virtualDevice virtualDevice to use.

Returns

Error Status. 0 on success.

Sets the number of analog channels in the datastream.

Parameters

NumberOfChannels	Number of analog channels.
------------------	----------------------------

Returns

Error Status. 0 on success.

Sets the number of analog channels in the datastream.

Parameters

NumberOfChannels	Number of analog channels.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

```
11.62.3.12 SetTriggerMaskValue() virtual uint32_t SetTriggerMaskValue (
    unsigned int mask,
    unsigned int value,
    unsigned int virtualDevice ) [virtual]
```

Defines a pattern on the digital dataword which will start a trigger when found.

Parameters

mask	Bits in the digital dataword which are monitored for a match with value.
------	--

value	Pattern which must match for the trigger to start.
-------	--

Returns

Error Status. 0 on success.

Sets the maximum number of samples per trigger.

Parameters

Returns

Error Status. 0 on success.

11.62.4 Property Documentation

11.62.4.1 AnalogGain virtual int AnalogGain [get]

The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

```
11.62.4.2 Gain virtual int Gain [get]
```

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.62.4.3 MeaAudioFunctionNet CMeaAudioFunctionNet^ MeaAudioFunctionNet [get]

11.62.4.4 MeaDigitalDataFunctionNet CMeaDigitalDataFunctionNet^ MeaDigitalDataFunctionNet [get]

```
11.62.4.5 MeaFeedbackFunctionNet CMeaFeedbackFunctionNet^ MeaFeedbackFunctionNet [get]
```


```
11.62.4.6 MeFunctionNet CMeFunctionNet^ MeFunctionNet [get]
```

```
11.62.4.7 W2100_FunctionNet CW2100_FunctionNet^ W2100_FunctionNet [get]
```

11.62.4.8 WClassicFunctionNet CWClassicFunctionNet^ WClassicFunctionNet [get]

11.63 CMeaDigitalDataFunctionNet Class Reference

Inheritance diagram for CMeaDigitalDataFunctionNet:

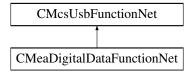

Public Member Functions

- CMeaDigitalDataFunctionNet (CMcsUsbNet[^] mcsusb)
- void SetDigitalData (unsigned int digital_value, unsigned int digital_value_mask)

Generate a value on the digital output.

• void SetDigitalData (unsigned int bit_number, bool value)

Generate a value on the digital output.

unsigned int GetDigitalData ()

Get the value of the digital output.

Additional Inherited Members

11.63.1 Constructor & Destructor Documentation

```
11.63.1.1 CMeaDigitalDataFunctionNet() [1/2] CMeaDigitalDataFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ meaDigitalFunctionPointerContainer)
```

11.63.1.2 CMeaDigitalDataFunctionNet() [2/2] CMeaDigitalDataFunctionNet (CMcsUsbNet^ mcsusb)

11.63.2 Member Function Documentation

11.63.2.1 GetDigitalData() unsigned int GetDigitalData ()

Get the value of the digital output.

Returns

Value on the digital data register.

Generate a value on the digital output.

Parameters

bit_number Bit number to change.

Parameters

value Bit value.

11.63.2.3 SetDigitalData() [2/2] void SetDigitalData (unsigned int digital_value, unsigned int digital_value_mask)

Generate a value on the digital output.

Parameters


digital_value Value to set.

Parameters

digital_value_mask	Mask for change.
--------------------	------------------

11.64 CMeaFeedbackFunctionNet Class Reference

Inheritance diagram for CMeaFeedbackFunctionNet:

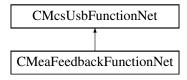

Public Member Functions

- CMeaFeedbackFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] meaFeedback
 FunctionNet)
- CMeaFeedbackFunctionNet (CMcsUsbNet[^] mcsusb)
- void FeedbackSetFeedback (unsigned char on, unsigned short digoutmask, unsigned short diginmask)
- unsigned int FeedbackGetSampleTimerCount ([System::Runtime::InteropServices::Out]unsigned int% CurrentCount, [System::Runtime::InteropServices::Out]unsigned int% LastKnownCount, [System::Runtime-::InteropServices::Out]bool% On)
- void FeedbackSetDigitalMapping (unsigned short channel, unsigned short outmapping, unsigned short inmapping)
- void FeedbackSetFilterParameter (unsigned char filter, array< short >[∧] parameters)
- void FeedbackSetFilterParameter32 (unsigned char filter, array< int >^ parameters)
- void FeedbackSetIIRFilterParameter (unsigned char filter, int length, array< double >^ parameters)
- void FeedbackSetMkFilter (unsigned char filter, String[^] filtertype, double cheb_ribble, String[^] passtype, int order, double alpha1, double alpha2)
- void FeedbackSetChannelFilter (short channel, char filter)
- void FeedbackSetGlobalChannelFilter (char filter, unsigned short firstchannel, unsigned short lastchannel)
- void FeedbackSetFilterOff ()
- void FeedbackSetNumberOfSpikeDetectors (unsigned short number)
- void FeedbackSetSpikeDetectorThreshold (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, int threshold1, int threshold2, short slope)
- void FeedbackSetNumberOfRateCounter (unsigned short number)
- void FeedbackSetRateCounter (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel)
- void FeedbackSetNumberOfRateDetectors (unsigned short number)
- void FeedbackSetRateDetector (unsigned short position, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, unsigned short pulses, unsigned int duration1, unsigned int duration2)
- void FeedbackSetNumberOfLogics (unsigned short number)
- void FeedbackSetLogic (unsigned short position, array< unsigned short >^ sourcechannel, unsigned short resultchannel, unsigned int lookup)
- void FeedbackSetNumberOfTriggers (unsigned short number)
- void FeedbackSetTrigger (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit)
- void FeedbackSetAnalogSource (AnalogSourceEnumNet AnalogSource, unsigned int Channels, unsigned int Offset)

Additional Inherited Members

11.64.1 Constructor & Destructor Documentation

```
CMcsUsbNet^ mcsusb,
            CMcsUsbFunctionPointerContainer<sup>∧</sup> meaFeedbackFunctionNet )
11.64.1.2 CMeaFeedbackFunctionNet() [2/2] CMeaFeedbackFunctionNet (
            CMcsUsbNet^ mcsusb )
11.64.2 Member Function Documentation
11.64.2.1 FeedbackGetSampleTimerCount() unsigned int FeedbackGetSampleTimerCount (
            [System::Runtime::InteropServices::Out] unsigned int% CurrentCount,
            [System::Runtime::InteropServices::Out] unsigned int% LastKnownCount,
            [System::Runtime::InteropServices::Out] bool% On )
11.64.2.2 FeedbackSetAnalogSource() void FeedbackSetAnalogSource (
            AnalogSourceEnumNet AnalogSource,
            unsigned int Channels,
            unsigned int Offset )
11.64.2.3 FeedbackSetChannelFilter() void FeedbackSetChannelFilter (
            short channel,
            char filter )
11.64.2.4 FeedbackSetDigitalMapping() void FeedbackSetDigitalMapping (
            unsigned short channel,
            unsigned short outmapping,
            unsigned short inmapping )
```

```
11.64.2.5 FeedbackSetFeedback() void FeedbackSetFeedback (
             unsigned char on,
             unsigned short digoutmask,
             unsigned short diginmask )
11.64.2.6 FeedbackSetFilterOff() void FeedbackSetFilterOff ()
11.64.2.7 FeedbackSetFilterParameter() void FeedbackSetFilterParameter (
             unsigned char filter,
             array < short >^{\wedge} parameters)
11.64.2.8 FeedbackSetFilterParameter32() void FeedbackSetFilterParameter32 (
             unsigned char filter,
             array< int >^{\wedge} parameters )
11.64.2.9 FeedbackSetGlobalChannelFilter() void FeedbackSetGlobalChannelFilter (
             char filter,
             unsigned short firstchannel,
             unsigned short lastchannel )
11.64.2.10 FeedbackSetIIRFilterParameter() void FeedbackSetIIRFilterParameter (
             unsigned char filter,
             int length,
             array< double >^{\land} parameters )
11.64.2.11 FeedbackSetLogic() void FeedbackSetLogic (
             unsigned short position,
             array< unsigned short >^{\wedge} sourcechannel,
             unsigned short resultchannel,
             unsigned int lookup )
```

```
11.64.2.12 FeedbackSetMkFilter() void FeedbackSetMkFilter (
              unsigned char filter,
              String^ filtertype,
              double cheb_ribble,
              String^{\wedge} passtype,
              int order,
              double alpha1,
              double alpha2 )
11.64.2.13 FeedbackSetNumberOfLogics() void FeedbackSetNumberOfLogics (
              unsigned short number )
11.64.2.14 FeedbackSetNumberOfRateCounter() void FeedbackSetNumberOfRateCounter (
              unsigned short number )
\textbf{11.64.2.15} \quad \textbf{FeedbackSetNumberOfRateDetectors()} \quad \texttt{void} \ \ \textbf{FeedbackSetNumberOfRateDetectors} \ \ \textbf{(}
              unsigned short number )
\textbf{11.64.2.16} \quad \textbf{FeedbackSetNumberOfSpikeDetectors()} \quad \texttt{void FeedbackSetNumberOfSpikeDetectors} \quad \textbf{(}
              unsigned short number )
11.64.2.17 FeedbackSetNumberOfTriggers() void FeedbackSetNumberOfTriggers (
              unsigned short number )
11.64.2.18 FeedbackSetRateCounter() void FeedbackSetRateCounter (
              unsigned short position,
              unsigned short sourcechannel,
              unsigned short resultchannel )
11.64.2.19 FeedbackSetRateDetector() void FeedbackSetRateDetector (
              unsigned short position,
              unsigned short resultchannel,
              unsigned short trigger,
              unsigned short totzeit,
              unsigned short pulses,
              unsigned int duration1,
              unsigned int duration2 )
```

11.64.2.20 FeedbackSetSpikeDetectorThreshold() void FeedbackSetSpikeDetectorThreshold (


```
unsigned short position,
unsigned short sourcechannel,
unsigned short resultchannel,
unsigned short trigger,
unsigned short totzeit,
int threshold1,
int threshold2,
short slope )
```

11.64.2.21 FeedbackSetTrigger() void FeedbackSetTrigger (

```
unsigned short position,
unsigned short sourcechannel,
unsigned short resultchannel,
unsigned short trigger,
unsigned short totzeit)
```

11.65 CMealmpedanceDeviceNet Class Reference

Inheritance diagram for CMealmpedanceDeviceNet:

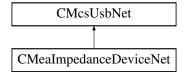

Public Member Functions

- CMealmpedanceDeviceNet ()
- CMeaImpedanceDeviceNet ()
- · void StartMeasurement (unsigned short channel)
- unsigned short GetReady ()
- unsigned short GetArraySize ()
- array< unsigned short > ^ GetResult ()
- unsigned short GetAdapterCode ()
- unsigned int GetImpedanceTestFrequency ()
- void SetImpedanceTestFrequency (unsigned int TestFrequency_Hertz)

Additional Inherited Members


11.65.1 Constructor & Destructor Documentation


```
11.65.1.2 ~CMealmpedanceDeviceNet() ~CMealmpedanceDeviceNet ()
11.65.2 Member Function Documentation
11.65.2.1 GetAdapterCode() unsigned short GetAdapterCode ( )
11.65.2.2 GetArraySize() unsigned short GetArraySize ()
11.65.2.3 GetImpedanceTestFrequency() unsigned int GetImpedanceTestFrequency ( )
11.65.2.4 GetReady() unsigned short GetReady ( )
11.65.2.5 GetResult() array<unsigned short> ^ GetResult ( )
11.65.2.6 SetImpedanceTestFrequency() void SetImpedanceTestFrequency (
            unsigned int TestFrequency_Hertz )
11.65.2.7 StartMeasurement() void StartMeasurement (
            unsigned short channel )
```

11.66 CMeasureTableDeviceNet Class Reference

CMeasureTableDeviceNet is the to control the MCS HLA device

 $Inheritance\ diagram\ for\ CMeasure Table Device Net:$

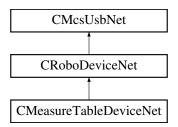

Public Member Functions

CMeasureTableDeviceNet (void)

Properties

• CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.66.1 Detailed Description

CMeasureTableDeviceNet is the to control the MCS HLA device

11.66.2 Constructor & Destructor Documentation

```
11.66.2.1 CMeasureTableDeviceNet() CMeasureTableDeviceNet ( void )
```


11.66.3 Property Documentation

```
11.66.3.1 Sensor CMcsBus_SensorNet^ Sensor [get]
```

11.67 CMeaSwitchDeviceNet Class Reference

The class to control the USB-MEA-Switch.

Inheritance diagram for CMeaSwitchDeviceNet:

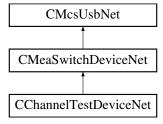

Public Member Functions

CMeaSwitchDeviceNet ()

Constructor.

∼CMeaSwitchDeviceNet ()

Destructor.

unsigned short GetNumber ()

Gets the number of boards in the device.

array< unsigned char > ^ GetPattern ()

Gets the pattern of the switches that are currently set in the device as char array.

array< bool > ^ GetPatternBool ()

Gets the pattern of the switches that are currently set in he device as bools.

void SetPattern (array< unsigned char >^ pattern)

Sets the pattern of switches from a char array.

void SetPatternBool (array< bool >^ pattern)

Sets the pattern of switches from a.

Additional Inherited Members

11.67.1 Detailed Description

The class to control the USB-MEA-Switch.

This class controls the settings of the USB-MEA-Switch. The box has two inputs for signals from a MEA amplifier. Each of the 64 outputs can be connected to one of the MEAs at the same channel.

11.67.2 Constructor & Destructor Documentation

```
11.67.2.1 CMeaSwitchDeviceNet() CMeaSwitchDeviceNet ()
```

Constructor.

$\textbf{11.67.2.2} \quad \sim \textbf{CMeaSwitchDeviceNet()} \quad \sim \texttt{CMeaSwitchDeviceNet()} \quad \\$

Destructor.

11.67.3 Member Function Documentation

```
11.67.3.1 GetNumber() unsigned short GetNumber ( )
```

Gets the number of boards in the device.

The MEA-Switch are delivered with 64 or 128 channels

```
11.67.3.2 GetPattern() array<unsigned char> ^{\wedge} GetPattern ( )
```

Gets the pattern of the switches that are currently set in the device as char array.

```
11.67.3.3 GetPatternBool() array<bool> ^ GetPatternBool ( )
```

Gets the pattern of the switches that are currently set in he device as bools.

```
11.67.3.4 SetPattern() void SetPattern ( array < unsigned char >^{\wedge} pattern )
```


Sets the pattern of switches from a char array.

Sets the pattern of switches from a.

11.68 CMeaUSBDeviceNet Class Reference

Class for data acquisition via ME and MEA USB amplifiers

Inheritance diagram for CMeaUSBDeviceNet:

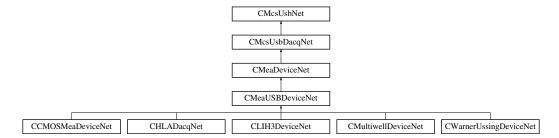

Public Member Functions

- CMeaUSBDeviceNet (OnChannelData[^] channelData, OnError[^] error)
 - Initializes a new instance of CMeaDeviceNet class.
- CMeaUSBDeviceNet ()

Initializes a new instance of CMeaDeviceNet class.

∼CMeaUSBDeviceNet ()

Additional Inherited Members

11.68.1 Detailed Description

Class for data acquisition via ME and MEA USB amplifiers

11.68.2 Constructor & Destructor Documentation

```
11.68.2.1 CMeaUSBDeviceNet() [1/2] CMeaUSBDeviceNet (
OnChannelData^ channelData,
OnError^ error )
```

Initializes a new instance of CMeaDeviceNet class.

Parameters

channelData	Handler to call when new data is available.
-------------	---

Parameters

error Handler to call when an error occurs.


11.68.2.2 CMeaUSBDeviceNet() [2/2] CMeaUSBDeviceNet ()

Initializes a new instance of CMeaDeviceNet class.

11.68.2.3 ~CMeaUSBDeviceNet() ~CMeaUSBDeviceNet ()

11.69 CMeFunctionNet Class Reference

Inheritance diagram for CMeFunctionNet:

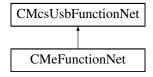

Public Member Functions

CMeFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] meFunctionPointer
 — Container)

Initializes a new instance of the CDacCalibrationFunctionNet class.

- CMeFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CMeFunctionNet (void)
- !CMeFunctionNet (void)
- · void SetTransformer (unsigned int index, bool onoff)

Additional Inherited Members

11.69.1 Detailed Description

11.69.2 Constructor & Destructor Documentation

Initializes a new instance of the CDacCalibrationFunctionNet class.

```
11.69.2.2 CMeFunctionNet() [2/2] CMeFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.69.2.3 ~CMeFunctionNet() virtual ~CMeFunctionNet (
void ) [virtual]
```


```
11.69.2.4 "!CMeFunctionNet() !CMeFunctionNet (
```

11.69.3 Member Function Documentation

11.70 CMultiBatteryChargerDeviceNet Class Reference

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

Inheritance diagram for CMultiBatteryChargerDeviceNet:

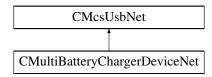

Public Member Functions

CMultiBatteryChargerDeviceNet ()

Initializes a new instance of the CMultiBatteryChargerDeviceNet class.

- virtual ~CMultiBatteryChargerDeviceNet ()
- !CMultiBatteryChargerDeviceNet ()
- uint32_t GetChargeCurrent (uint32_t NrChannel)

gets the charge current; unit: mA

uint32_t GetDischargeCurrent (uint32_t NrChannel)

gets the discharge current; unit: mA

void SetDischargeCurrentSetPoint (uint32_t NrChannel, uint32_t DischargeCurrent_mA)

sets the setpoint for the discharge current; unit: mA

uint32_t GetDischargeCurrentSetPoint (uint32_t NrChannel)

gets the setpoint for the discharge current; unit: mA

• void SetFinalDischargeVoltage (uint32_t NrChannel, uint32_t FinalDischargeVoltage_mV)

sets the final discharge voltage; unit: mV

uint32_t GetFinalDischargeVoltage (uint32_t NrChannel)

gets the final discharge voltage; unit: mV

uint32_t GetDischargeCapacity (uint32_t NrChannel)

gets the discharge capacity; unit: ?Ah

uint32_t GetChargeCapacity (uint32_t NrChannel)

gets the charge capacity; unit: ?Ah

uint32_t GetBatteryVoltage (uint32_t NrChannel)

gets the battery voltage; unit: mV

uint32_t GetChannels ()

gets number of channels

void SetRatedCapacityVolatile (uint32_t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)

sets the rated capacity (i.e. charge current) without storing it persistently

void SetChargingMode (uint32 t NrChannel, MbcChargingModeEnumNet NewOperatingMode)

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

• MbcChargingModeEnumNet GetChargingMode (uint32_t NrChannel)

gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

MbcChannelStateEnumNet GetChannelState (uint32_t NrChannel)

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

void CapacityTest (uint32_t NrChannel)

start capacity test on channel

· void ChannelReset (uint32 t NrChannel)

cancel charging and capacity test functions; check if battery is connected

- void SetChargingPCoefficient (uint32_t pCoefficient)
 sets the p-coefficient for charging in mA/V / nominal charging current
- uint32 t GetChargingPCoefficient ()

gets the p-coefficient for charging in mA/V / nominal charging current

void SetRatedCapacity (uint32_t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)
 sets the rated capacity

 MbcRatedCapacityEnumNet GetRatedCapacity (uint32_t NrChannel) gets the rated capacity

Additional Inherited Members

11.70.1 Detailed Description

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

11.70.2 Constructor & Destructor Documentation

11.70.2.1 CMultiBatteryChargerDeviceNet() CMultiBatteryChargerDeviceNet ()

Initializes a new instance of the CMultiBatteryChargerDeviceNet class.

```
11.70.2.2 ~CMultiBatteryChargerDeviceNet() virtual ~CMultiBatteryChargerDeviceNet ( ) [virtual]
```

```
11.70.2.3 "!CMultiBatteryChargerDeviceNet() !CMultiBatteryChargerDeviceNet ( )
```

11.70.3 Member Function Documentation

start capacity test on channel

Parameters

NrChannel | the channel number

```
11.70.3.2 ChannelReset() void ChannelReset ( uint32_t NrChannel)
```

cancel charging and capacity test functions; check if battery is connected

Parameters

NrChannel the channel number

11.70.3.3 GetBatteryVoltage() uint32_t GetBatteryVoltage (uint32_t NrChannel)

gets the battery voltage; unit: mV

Parameters

Returns

the battery voltage in mV

11.70.3.4 GetChannels() uint32_t GetChannels ()

gets number of channels

Returns

number of channels

11.70.3.5 **GetChannelState()** MbcChannelStateEnumNet GetChannelState (uint32_t NrChannel)

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

NrChannal	the channel number
NrChannel	the channel number

Returns

the current state

```
11.70.3.6 GetChargeCapacity() uint32_t GetChargeCapacity ( uint32_t NrChannel)
```

gets the charge capacity; unit: ?Ah

Parameters

NrChannel the channel number	
------------------------------	--

Returns

the capacity in uAh

gets the charge current; unit: mA

Parameters

NrChannel	the channel number
i ii Onamo	tilo olialilloi lialillooi

Returns

the measured charge current in mA

11.70.3.8 GetChargingMode() MbcChargingModeEnumNet GetChargingMode (uint32_t NrChannel)

gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

NrChannel the channel number

Returns

the charging mode

11.70.3.9 GetChargingPCoefficient() uint32_t GetChargingPCoefficient ()

gets the p-coefficient for charging in mA/V / nominal charging current

Returns

the p-coefficient

11.70.3.10 GetDischargeCapacity() uint32_t GetDischargeCapacity (uint32_t NrChannel)

gets the discharge capacity; unit: ?Ah

Parameters

NrChannel the channel number

Returns

the capacity in uAh

11.70.3.11 GetDischargeCurrent() uint32_t GetDischargeCurrent (uint32_t NrChannel)

gets the discharge current; unit: mA

Parameters

NrChannel | the channel number

Returns

the measured discharge current in mA

11.70.3.12 GetDischargeCurrentSetPoint() uint32_t GetDischargeCurrentSetPoint (uint32_t NrChannel)

gets the setpoint for the discharge current; unit: mA

Parameters

NrChannel the channel number

Returns

the discharge current in mA

```
11.70.3.13 GetFinalDischargeVoltage() uint32_t GetFinalDischargeVoltage ( uint32_t NrChannel )
```

gets the final discharge voltage; unit: mV

Parameters

NrChannel the channel number	oer
------------------------------	-----

Returns

the battery voltage in mV at the end of discharge

11.70.3.14 **GetRatedCapacity()** MbcRatedCapacityEnumNet GetRatedCapacity (uint32_t NrChannel)

gets the rated capacity

Parameters

NrChannel	the channel number
i ii Onamo	tilo olialilloi lialillooi

Returns

the capacity

```
11.70.3.15 SetChargingMode() void SetChargingMode (
```

uint32_t NrChannel,
MbcChargingModeEnumNet NewOperatingMode)

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

NrChannel	the channel number
NewOperatingMode	the charging mode

11.70.3.16 SetChargingPCoefficient() void SetChargingPCoefficient (uint32_t pCoefficient)

sets the p-coefficient for charging in mA/V / nominal charging current

Parameters

pCoefficient	the p-coefficient
--------------	-------------------

sets the setpoint for the discharge current; unit: mA

Parameters

NrChannel	the channel number
DischargeCurrent_mA	the discharge current in mA

sets the final discharge voltage; unit: mV

Parameters

NrChannel	the channel number
FinalDischargeVoltage_mV	the battery voltage in mV at the end of discharge

11.70.3.19 SetRatedCapacity() void SetRatedCapacity (uint32_t NrChannel,

 ${\tt MbcRatedCapacityEnumNet}\ \textit{NewRatedCapacity}\)$

sets the rated capacity

NrChannel	the channel number
NewRatedCapacity	the capacity

sets the rated capacity (i.e. charge current) without storing it persistently


Parameters

NrChannel	the channel number
NewRatedCapacity	the capacity

11.71 CMultiwellCallbackFunctionNet Class Reference

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

Inheritance diagram for CMultiwellCallbackFunctionNet:

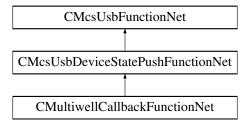

Public Member Functions

- delegate void OnGetPlateClampStateByHeadstage (uint32_t Headstage, PlateClampEnumNet plateState)
- CMultiwellCallbackFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pMultiwell←
 CallbackFunctionPointerContainer)

Initializes a new instance of the CMultiwellCallbackFunctionNet class.

- CMultiwellCallbackFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CMultiwellCallbackFunctionNet ()
- !CMultiwellCallbackFunctionNet ()
- PlateClampEnumNet GetPlateClampStateByHeadstage (uint32_t Headstage)

Gets the state of the plate

Events

Event fires when the plate state for the headstage number has changed

Additional Inherited Members

11.71.1 Detailed Description

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

11.71.2 Constructor & Destructor Documentation

```
11.71.2.1 CMultiwellCallbackFunctionNet() [1/2] CMultiwellCallbackFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pMultiwellCallbackFunctionPointerContainer)
```

Initializes a new instance of the CMultiwellCallbackFunctionNet class.

```
11.71.2.2 CMultiwellCallbackFunctionNet() [2/2] CMultiwellCallbackFunctionNet ( CMcsUsbNet^ mcsusb )
```

```
11.71.2.3 ~CMultiwellCallbackFunctionNet() virtual ~CMultiwellCallbackFunctionNet ( ) [virtual]
```

```
11.71.2.4 "!CMultiwellCallbackFunctionNet() !CMultiwellCallbackFunctionNet ()
```

11.71.3 Member Function Documentation

11.71.3.1 GetPlateClampStateByHeadstage() PlateClampEnumNet GetPlateClampStateByHeadstage (uint32_t Headstage)

Gets the state of the plate

Parameters

Headstage	The headstage number
-----------	----------------------

Returns

The plate state

11.71.4 Event Documentation


11.71.4.1 GetPlateClampStateByHeadstageEvent OnGetPlateClampStateByHeadstage[∧] GetPlateClamp↔ StateByHeadstageEvent [add], [remove], [raise]

Event fires when the plate state for the headstage number has changed

11.72 CMultiwellDeviceNet Class Reference

CMultiwellDeviceNet is the class to access the Multiwell device.

Inheritance diagram for CMultiwellDeviceNet:

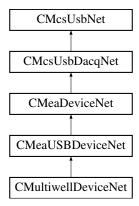

Public Member Functions

• CMultiwellDeviceNet ()

Initializes a new instance of the CMultiwellDeviceNet class.

- virtual ∼CMultiwellDeviceNet ()
- !CMultiwellDeviceNet ()
- PlateClampEnumNet GetPlateClampState ()

Gets the state of the Multiwell plate clamp.

• PlateClampEnumNet GetPlateClampState (uint32_t Headstage)

Gets the state of the plate

void OpenPlateClamp ()

Opens the plate clamp.

void ClosePlateClamp ()

Closes the plate clamp.

void StopPlateClamp ()

Stops the plate clamp movement.

uint32_t GetPlateClampLockState ()

Gets the state of the plate clamp lock.

void LockPlateClamp ()

Locks the plate clamp.

void UnlockPlateClamp ()

Unlocks the plate clamp.

• MultiwellPlateTypeEnumNet GetPlateType ()

Gets the plate type.

MultiwellPlateTypeEnumNet GetPlateType (uint32_t Headstage)

Gets the plate type.

void SetPlateType (MultiwellPlateTypeEnumNet plateType)

Sets the plate type.

• void SetPlateType (uint32_t Headstage, MultiwellPlateTypeEnumNet plateType)

Sets the plate type.

void SetPlateMux (uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

void SetPlateMux (uint32_t Headstage, uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

uint32 t GetPlateMux ()

Gets the selected quarter of the electrodes on a high density Multiwell plate.

uint32 t GetPlateMux (uint32 t Headstage)

Gets the selected quarter of the electrodes on a high density Multiwell plate.

bool IsPlateTypeValid ()

Checks whether the plate type is valid, meaning all pins have contact.

bool IsPlateTypeValid (uint32_t Headstage)

Checks whether the plate type is valid, meaning all pins have contact.

void SetPowerMuxPlate (uint32_t Headstage, bool powerOn)

On the Multiwell Mini device, turn Power to the MUX Plate On or Off.

bool GetPowerMuxPlate (uint32_t Headstage)

On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.

Additional Inherited Members

11.72.1 Detailed Description

CMultiwellDeviceNet is the class to access the Multiwell device.

11.72.2 Constructor & Destructor Documentation

11.72.2.1 CMultiwellDeviceNet() CMultiwellDeviceNet ()

Initializes a new instance of the CMultiwellDeviceNet class.

$\textbf{11.72.2.2} \quad \sim \textbf{CMultiwellDeviceNet()} \quad \text{virtual} \quad \sim \texttt{CMultiwellDeviceNet ()} \quad [\texttt{virtual}]$

11.72.2.3 "!CMultiwellDeviceNet() !CMultiwellDeviceNet ()

11.72.3 Member Function Documentation

11.72.3.1 ClosePlateClamp() void ClosePlateClamp ()

Closes the plate clamp.

11.72.3.2 GetPlateClampLockState() uint32_t GetPlateClampLockState ()

Gets the state of the plate clamp lock.

Returns

the state of the plate lock (unlocked/locked)

11.72.3.3 GetPlateClampState() [1/2] PlateClampEnumNet GetPlateClampState ()

Gets the state of the Multiwell plate clamp.

Returns

the state of the plate clamp (open/closed)

11.72.3.4 GetPlateClampState() [2/2] PlateClampEnumNet GetPlateClampState (uint32_t Headstage)

Gets the state of the plate

Parameters

Headstage	The headstage number

Returns

The plate state

11.72.3.5 GetPlateMux() [1/2] uint32_t GetPlateMux ()

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Returns

the selected quarter

```
11.72.3.6 GetPlateMux() [2/2] uint32_t GetPlateMux ( uint32_t Headstage )
```

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

the selected quarter

```
11.72.3.7 GetPlateType() [1/2] MultiwellPlateTypeEnumNet GetPlateType ( )
```

Gets the plate type.

Returns

the plate type

```
11.72.3.8 GetPlateType() [2/2] MultiwellPlateTypeEnumNet GetPlateType ( uint32_t Headstage )
```

Gets the plate type.

Parameters

```
Headstage The headstage to query.
```

Returns

the plate type

```
11.72.3.9 GetPowerMuxPlate() bool GetPowerMuxPlate ( uint32_t Headstage )
```

On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.

Parameters

Headstage The headstage to que

Returns

"true" Power is On, "false" Power is Off

11.72.3.10 IsPlateTypeValid() [1/2] bool IsPlateTypeValid ()

Checks whether the plate type is valid, meaning all pins have contact.

Returns

"true" when all pins have contact, otherwise "false".

11.72.3.11 **IsPlateTypeValid()** [2/2] bool IsPlateTypeValid (uint32_t *Headstage*)

Checks whether the plate type is valid, meaning all pins have contact.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

"true" when all pins have contact, otherwise "false".

11.72.3.12 LockPlateClamp() void LockPlateClamp ()

Locks the plate clamp.

11.72.3.13 OpenPlateClamp() void OpenPlateClamp ()

Opens the plate clamp.

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

Headstage	The headstage to query.
muxSelection	the selected quarter

11.72.3.15 SetPlateMux() [2/2] void SetPlateMux (uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

muxSelection	the selected quarter
--------------	----------------------

Sets the plate type.

Parameters

```
plateType the plate type
```

Sets the plate type.

Parameters

Headstage	The headstage to query.
plateType	the plate type

On the Multiwell Mini device, turn Power to the MUX Plate On or Off.

Headstage	The headstage to query.	
powerOn	"true" to turn Power On, "false" to turn Power Off	

11.72.3.19 StopPlateClamp() void StopPlateClamp ()

Stops the plate clamp movement.


11.72.3.20 UnlockPlateClamp() void UnlockPlateClamp ()

Unlocks the plate clamp.

11.73 CMultiwellOptoStimFunctionNet Class Reference

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

Inheritance diagram for CMultiwellOptoStimFunctionNet:

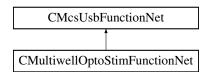

Public Member Functions

CMultiwellOptoStimFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pMultiwell←
 OptoStimFunctionPointerContainer)

Initializes a new instance of the CMultiwellOptoStimFunctionNet class.

- CMultiwellOptoStimFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CMultiwellOptoStimFunctionNet ()
- !CMultiwellOptoStimFunctionNet ()
- uint32_t GetWaveLengthInNanometer (uint16_t channel)
- uint32_t GetAbsMaxCurrentInMicroAmp (uint16_t channel)
- uint32_t GetMaxDurationHighCurrentInMicroSec (uint16_t channel)
- uint32_t GetMaxDutyCycleHighCurrent (uint16_t channel)
- uint32_t GetPermanentCurrentInMicroAmp (uint16_t channel)
- uint32_t GetColorRgb (uint16_t channel)
- String ^ GetColorStr (uint16_t channel)
- void SetWaveLengthInNanometer (uint16_t channel, uint32_t WaveLength_nm)
- void SetAbsMaxCurrentInMicroAmp (uint16_t channel, uint32_t AbsoluteMaxCurrent_uA)
- void SetMaxDurationHighCurrentInMicroSec (uint16_t channel, uint32_t AbsoluteMaxDuration_us)
- void SetMaxDutyCycleHighCurrent (uint16 t channel, uint32 t MaxDutyCycleHighCurrent)
- void SetPermanentCurrentInMicroAmp (uint16_t channel, uint32_t PermanentCurrent_uA)
- void SetColorRgb (uint16_t channel, uint32_t ColorRGB)
- void SetColorStr (uint16_t channel, String[^] ColorString)

Additional Inherited Members

11.73.1 Detailed Description

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

11.73.2 Constructor & Destructor Documentation

```
11.73.2.1 CMultiwellOptoStimFunctionNet() [1/2] CMultiwellOptoStimFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pMultiwellOptoStimFunctionPointerContainer )
```

Initializes a new instance of the CMultiwellOptoStimFunctionNet class.

```
11.73.2.2 CMultiwellOptoStimFunctionNet() [2/2] CMultiwellOptoStimFunctionNet (
CMcsUsbNet^ mcsusb)
```

11.73.2.3 ~CMultiwellOptoStimFunctionNet() virtual ~CMultiwellOptoStimFunctionNet () [virtual]

```
11.73.2.4 "!CMultiwellOptoStimFunctionNet() !CMultiwellOptoStimFunctionNet ( )
```

11.73.3 Member Function Documentation

```
11.73.3.1 GetAbsMaxCurrentInMicroAmp() uint32_t GetAbsMaxCurrentInMicroAmp ( uint16_t channel)
```

Parameters

```
channel the (analog) channel number
```

Returns

absolute max. current; unit: uA

```
11.73.3.2 GetColorRgb() uint32_t GetColorRgb ( uint16_t channel)
```

channel	the (analog) channel number
---------	-----------------------------

Returns

RGB-value of LED color

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

LED color as string

11.73.3.4 **GetMaxDurationHighCurrentInMicroSec()** uint32_t GetMaxDurationHighCurrentInMicroSec (uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

max. duration the LED can stand the abs. max current; unit: us

11.73.3.5 GetMaxDutyCycleHighCurrent() uint32_t GetMaxDutyCycleHighCurrent (uint16_t *channel*)

Parameters

channel	the (analog) channel number

max. duty cycle at max. current; unit: 100*%

11.73.3.6 GetPermanentCurrentInMicroAmp() uint32_t GetPermanentCurrentInMicroAmp (

uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

max. current the LED can stand when always switched on; unit: uA

11.73.3.7 GetWaveLengthInNanometer() uint32_t GetWaveLengthInNanometer (

uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

wavelength of this channel's LEDs; unit: nm

$\textbf{11.73.3.8} \quad \textbf{SetAbsMaxCurrentInMicroAmp()} \quad \texttt{void SetAbsMaxCurrentInMicroAmp} \quad \textbf{(}$

```
uint16_t channel,
uint32_t AbsoluteMaxCurrent_uA )
```

Parameters

channel	the (analog) channel number
AbsoluteMaxCurrent_uA	absolute max. current; unit: uA

11.73.3.9 SetColorRgb() void SetColorRgb (

```
uint16_t channel,
uint32_t ColorRGB )
```

channel	the (analog) channel number
ColorRGB	RGB-value of LED color

Parameters

channel	the (analog) channel number
ColorString	LED color as string

Parameters

channel	the (analog) channel number	
AbsoluteMaxDuration_us	max. duration the LED can stand the abs. max current; unit: us	

Parameters

channel	the (analog) channel number
MaxDutyCycleHighCurrent	max. duty cycle at max. current; unit: 100*%

Parameters


channel	the (analog) channel number	
PermanentCurrent_uA	max. current the LED can stand when always switched on; unit: uA	

Parameters

channel	the (analog) channel number
WaveLength_nm	wavelength of this channel's LEDs; unit: nm

11.74 CNF_GenDeviceNet Class Reference

Inheritance diagram for CNF_GenDeviceNet:

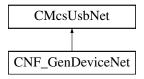

Public Member Functions

- CNF_GenDeviceNet (void)
- \sim CNF_GenDeviceNet (void)
- void Set_Values (unsigned int frequency, unsigned int amplitude)

Additional Inherited Members

11.74.1 Constructor & Destructor Documentation

```
11.74.1.1 CNF_GenDeviceNet() CNF_GenDeviceNet ( void )
```


11.74.1.2
$$\sim$$
 CNF_GenDeviceNet() \sim CNF_GenDeviceNet (void)

11.74.2 Member Function Documentation

```
11.74.2.1 Set_Values() void Set_Values (
          unsigned int frequency,
          unsigned int amplitude )
```

11.75 COctoPotDeviceNet Class Reference

Inheritance diagram for COctoPotDeviceNet:

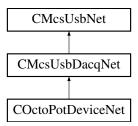

Public Member Functions

- COctoPotDeviceNet (void)
- COctoPotDeviceNet (OnChannelData^ channelData, OnError^ error)
- uint32_t SetOutputRate (uint32_t rate)
- uint32_t SetBathclamp (unsigned int block, bool enable)
- uint32 t SetDacValue (int channel, int value)
- uint32_t SetDacAutoControl (unsigned int channel)
- uint32_t SetPidParameter (unsigned int channel, int const_p, int const_i, int shift_p, int shift_i)
- uint32_t SetRampParameter (unsigned int channel, int start, int min, int max, int slope, int slope2, int pause, unsigned int samples)
- uint32_t RampStart (int channelmap)
- uint32_t SetSineParameter (unsigned int channel, int amplitude)
- uint32_t SineStart (int channelmap)
- uint32_t SetPatternListEntry (unsigned int channel, unsigned int position, unsigned int duration, int value)
- uint32_t PatternListStart (int channelmap)
- uint32_t SetAdcOffset (unsigned int channel, int offset)
- · uint32_t SetDacOffset (unsigned int channel, int offset)
- uint32_t ResetAdcOffset (unsigned int channel)
- uint32_t ResetDacOffset (unsigned int channel)
- uint32 t BurnAdcOffset ()
- uint32_t BurnDacOffset ()
- uint32_t GetAdcOffset (unsigned int channel, [System::Runtime::InteropServices::Out] int [∧] offset)
- uint32 t GetDacOffset (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- uint32 t SetAmplificationSwitch (unsigned int channel, unsigned int state)
- uint32_t SetChannelSwitch (unsigned int channel, unsigned int state)
- uint32_t SetNumberOfChannels (unsigned int NumberOfChannels)
- uint32 t EnableDigitalIn (bool enable)
- uint32_t EnableTimestamp (bool enable)
- uint32 t EnableChecksum (bool enable)

Additional Inherited Members

11.75.1 Constructor & Destructor Documentation


```
11.75.1.1 COctoPotDeviceNet() [1/2] COctoPotDeviceNet (
             void )
11.75.1.2 COctoPotDeviceNet() [2/2] COctoPotDeviceNet (
             OnChannelData^ channelData,
             OnError<sup>∧</sup> error )
11.75.2 Member Function Documentation
11.75.2.1 BurnAdcOffset() uint32_t BurnAdcOffset ( )
11.75.2.2 BurnDacOffset() uint32_t BurnDacOffset ( )
11.75.2.3 EnableChecksum() uint32_t EnableChecksum (
             bool enable )
11.75.2.4 EnableDigitalIn() uint32_t EnableDigitalIn (
             bool enable )
11.75.2.5 EnableTimestamp() uint32_t EnableTimestamp (
             bool enable )
11.75.2.6 GetAdcOffset() uint32_t GetAdcOffset (
             unsigned int channel,
             [System::Runtime::InteropServices::Out] int ^ offset )
11.75.2.7 GetDacOffset() uint32_t GetDacOffset (
             unsigned int channel,
             [System::Runtime::InteropServices::Out] int ^{\land} offset )
```

```
11.75.2.8 PatternListStart() uint32_t PatternListStart (
              int channelmap )
11.75.2.9 RampStart() uint32_t RampStart (
              int channelmap )
\textbf{11.75.2.10} \quad \textbf{ResetAdcOffset()} \quad \texttt{uint32\_t} \;\; \texttt{ResetAdcOffset} \;\; \texttt{(}
              unsigned int channel )
11.75.2.11 ResetDacOffset() uint32_t ResetDacOffset (
              unsigned int channel )
11.75.2.12 SetAdcOffset() uint32_t SetAdcOffset (
              unsigned int channel,
              int offset )
11.75.2.13 SetAmplificationSwitch() uint32_t SetAmplificationSwitch (
              unsigned int channel,
              unsigned int state )
11.75.2.14 SetBathclamp() uint32_t SetBathclamp (
              unsigned int block,
              bool enable )
11.75.2.15 SetChannelSwitch() uint32_t SetChannelSwitch (
              unsigned int channel,
              unsigned int state )
11.75.2.16 SetDacAutoControl() uint32_t SetDacAutoControl (
              unsigned int channel )
```

```
11.75.2.17 SetDacOffset() uint32_t SetDacOffset (
             unsigned int channel,
             int offset )
11.75.2.18 SetDacValue() uint32_t SetDacValue (
             int channel,
             int value )
11.75.2.19 SetNumberOfChannels() uint32_t SetNumberOfChannels (
             unsigned int NumberOfChannels )
11.75.2.20 SetOutputRate() uint32_t SetOutputRate (
             uint32_t rate )
11.75.2.21 SetPatternListEntry() uint32_t SetPatternListEntry (
             unsigned int channel,
             unsigned int position,
             unsigned int duration,
             int value )
11.75.2.22 SetPidParameter() uint32_t SetPidParameter (
             unsigned int channel,
             int const_p,
             int const_i,
             int shift_p,
             int shift_i )
11.75.2.23 SetRampParameter() uint32_t SetRampParameter (
             unsigned int channel,
             int start,
             int min,
             int max,
             int slope,
             int slope2,
             int pause,
             unsigned int samples )
```

11.76 COkuvisionStimulatorDeviceNet Class Reference

Inheritance diagram for COkuvisionStimulatorDeviceNet:

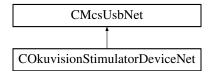

Public Member Functions

- · COkuvisionStimulatorDeviceNet (void)
- ~COkuvisionStimulatorDeviceNet (void)
- void SetPulseform (int channel, int current, int pulsewidth, int periode, int duration)
- void SetMaxPower (int channel, int power)
- int GetMaxPower (int channel)
- void SetMaxVoltage (int channel, int voltage)
- int GetMaxVoltage (int channel)
- void SetCheckVoltage (int channel, int voltage)
- int GetCheckVoltage (int channel)
- int GetVoltage (int channel)
- · void SetDACOffset (int channel, int part, int offset)
- int GetDACOffset (int channel, int part)
- void SetRTC (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)
- void GetRTC ([System::Runtime::InteropServices::Out] uint8_t% year, [System::Runtime::InteropServices
 ::Out] uint8_t% month, [System::Runtime::InteropServices::Out] uint8_t% day, [System::Runtime::Interop
 Services::Out] uint8_t% hour, [System::Runtime::InteropServices::Out] uint8_t% minute, [System::Runtime
 ::InteropServices::Out] uint8_t% second)
- void SetRTC (DateTime timestamp)
- DateTime GetRTC ()
- void GetStimulatorStatus ([System::Runtime::InteropServices::Out] int% startstop, [System::Runtime::

 InteropServices::Out] int% last_error, [System::Runtime::InteropServices::Out] int% battery_status)
- void SetCurrentFactor (int channel, int factor)
- · int GetCurrentFactor (int channel)

Additional Inherited Members

11.76.1 Constructor & Destructor Documentation

```
11.76.1.1 COkuvisionStimulatorDeviceNet() COkuvisionStimulatorDeviceNet (
            void )
11.76.1.2 ~COkuvisionStimulatorDeviceNet() ~COkuvisionStimulatorDeviceNet (
             void )
11.76.2 Member Function Documentation
11.76.2.1 GetCheckVoltage() int GetCheckVoltage (
             int channel )
11.76.2.2 GetCurrentFactor() int GetCurrentFactor (
             int channel )
11.76.2.3 GetDACOffset() int GetDACOffset (
            int channel,
             int part )
11.76.2.4 GetMaxPower() int GetMaxPower (
            int channel )
11.76.2.5 GetMaxVoltage() int GetMaxVoltage (
             int channel )
```


```
 \textbf{11.76.2.6} \quad \textbf{GetPulseform()} \quad \texttt{void GetPulseform ()} \\
              int channel,
              [System::Runtime::InteropServices::Out] int% current,
              [System::Runtime::InteropServices::Out] int% pulsewidth,
              [System::Runtime::InteropServices::Out] int% periode,
              [System::Runtime::InteropServices::Out] int% duration )
11.76.2.7 GetRTC() [1/2] DateTime GetRTC ( )
11.76.2.8 GetRTC() [2/2] void GetRTC (
              [System::Runtime::InteropServices::Out] uint8_t% year,
              [System::Runtime::InteropServices::Out] uint8_t% month,
              [System::Runtime::InteropServices::Out] uint8_t% day,
              [System::Runtime::InteropServices::Out] uint8_t% hour,
              [System::Runtime::InteropServices::Out] uint8_t% minute,
              [System::Runtime::InteropServices::Out] uint8_t% second )
11.76.2.9 GetStimulatorStatus() void GetStimulatorStatus (
              [System::Runtime::InteropServices::Out] int% startstop,
              [{\tt System::Runtime::InteropServices::Out}] \  \, {\tt int} \\ \textit{ last\_error,}
              [System::Runtime::InteropServices::Out] int% battery_status )
11.76.2.10 GetVoltage() int GetVoltage (
              int channel )
11.76.2.11 SetCheckVoltage() void SetCheckVoltage (
              int channel,
              int voltage )
11.76.2.12 SetCurrentFactor() void SetCurrentFactor (
              int channel,
              int factor )
```

```
11.76.2.13 SetDACOffset() void SetDACOffset (
             int channel,
             int part,
             int offset )
11.76.2.14 SetMaxPower() void SetMaxPower (
             int channel,
             int power )
11.76.2.15 SetMaxVoltage() void SetMaxVoltage (
             int channel,
             int voltage )
11.76.2.16 SetPulseform() void SetPulseform (
             int channel,
             int current,
             int pulsewidth,
             int periode,
             int duration )
11.76.2.17 SetRTC() [1/2] void SetRTC (
             DateTime timestamp )
11.76.2.18 SetRTC() [2/2] void SetRTC (
             uint8_t year,
             uint8_t month,
             uint8_t day,
             uint8_t hour,
             uint8_t minute,
             uint8_t second )
```

11.77 CPatchServerDeviceNet Class Reference

CPatchServerDeviceNet is the class to control the MCS PatchServer device

Inheritance diagram for CPatchServerDeviceNet:

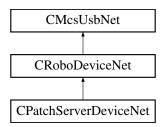

Public Member Functions

CPatchServerDeviceNet (void)

Properties

• CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.77.1 Detailed Description

CPatchServerDeviceNet is the class to control the MCS PatchServer device

11.77.2 Constructor & Destructor Documentation


```
11.77.2.1 CPatchServerDeviceNet() CPatchServerDeviceNet (
void )
```

11.77.3 Property Documentation

```
11.77.3.1 Sensor CMcsBus_SensorNet^ Sensor [get]
```

11.78 CPathIdentDeviceNet Class Reference

Inheritance diagram for CPathIdentDeviceNet:

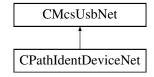

Public Member Functions

- CPathIdentDeviceNet (void)
- ∼CPathIdentDeviceNet (void)
- void Set_Values (unsigned int frequency, unsigned int amplitude)
- void Measure ([System::Runtime::InteropServices::Out] unsigned int% phase, [System::Runtime::Interop←
 Services::Out] unsigned int% amplitude)

Additional Inherited Members

11.78.1 Constructor & Destructor Documentation


```
11.78.1.1 CPathIdentDeviceNet() CPathIdentDeviceNet (
             void )
11.78.1.2 ~ CPathIdentDeviceNet() ~ CPathIdentDeviceNet (
             void )
11.78.2 Member Function Documentation
```

```
11.78.2.1 Measure() void Measure (
             [System::Runtime::InteropServices::Out] unsigned int% phase,
             [System::Runtime::InteropServices::Out] unsigned int% amplitude )
11.78.2.2 Set_Values() void Set_Values (
             unsigned int frequency,
```

11.79 CPedoterDeviceNet Class Reference

unsigned int amplitude)

Inheritance diagram for CPedoterDeviceNet:

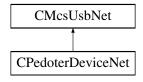

Public Member Functions

CPedoterDeviceNet ()

Initializes a new instance of the CPedoterDeviceNet class.

- virtual ∼CPedoterDeviceNet ()
- !CPedoterDeviceNet ()
- uint32_t GetCommand (uint16_t Argument)

Get value from the pedoter device

void SetCommand (uint16_t Argument, uint32_t pData)

Set value on the pedoter device

Additional Inherited Members

11.79.1 Detailed Description

11.79.2 Constructor & Destructor Documentation

```
11.79.2.1 CPedoterDeviceNet() CPedoterDeviceNet ()
```

Initializes a new instance of the CPedoterDeviceNet class.

```
11.79.2.2 ~ CPedoterDeviceNet() virtual ~ CPedoterDeviceNet () [virtual]
```

```
11.79.2.3 "!CPedoterDeviceNet() !CPedoterDeviceNet ()
```

11.79.3 Member Function Documentation

```
11.79.3.1 GetCommand() uint32_t GetCommand ( uint16_t Argument )
```

Get value from the pedoter device

Parameters

```
Argument argument
```

Returns

value


Set value on the pedoter device

Argument	argument	
pData	value	

11.80 CPeristalticPumpDeviceNet Class Reference

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

Inheritance diagram for CPeristalticPumpDeviceNet:

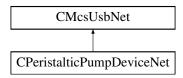

Public Member Functions

- CPeristalticPumpDeviceNet (void)

 Initialize a new instance of the CPeristalticPumpDeviceNet class.
- ~CPeristalticPumpDeviceNet (void)

Properties

CMcsBus_MotorControlNet[^] McsBus_MotorControl [get]

Additional Inherited Members

11.80.1 Detailed Description

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

11.80.2 Constructor & Destructor Documentation

11.80.2.1 CPeristalticPumpDeviceNet() CPeristalticPumpDeviceNet (void)

Initialize a new instance of the CPeristalticPumpDeviceNet class.


```
11.80.2.2 ~CPeristalticPumpDeviceNet() ~CPeristalticPumpDeviceNet (
```

11.80.3 Property Documentation

11.80.3.1 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.81 CPgaDeviceNet Class Reference

Inheritance diagram for CPgaDeviceNet:

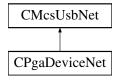

Public Member Functions

- · CPgaDeviceNet ()
- ∼CPgaDeviceNet ()
- uint32_t GetNumFrequencyRanges ([System::Runtime::InteropServices::Out]int% numRanges)
- uint32_t GetFrequencyRange (int rangeIndex, [System::Runtime::InteropServices::Out]int% low, [System::

 Runtime::InteropServices::Out]int% high, [System::Runtime::InteropServices::Out]int% channels, [System
 ::Runtime::InteropServices::Out]int% gain)
- uint32_t GetNumAmplifications ([System::Runtime::InteropServices::Out]int% number)
- uint32_t GetAmplification (int index, [System::Runtime::InteropServices::Out]int% amplification, [System:: \leftarrow Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)
- uint32_t DefineNumFrequencyRanges (int rnum)
- uint32_t DefineFrequencyRange (int index, int low, int high, int channels, int gain)
- uint32_t DefineNumAmplifications (int number)
- uint32_t DefineAmplification (int index, int amplification, int poti1, int poti2)
- uint32_t SetGain (int channel, int Gain, int poti1, int poti2)
- uint32_t GetGain (int channel, [System::Runtime::InteropServices::Out]int% Gain, [System::Runtime::
 InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)
- uint32 t ApplyGains ()

Additional Inherited Members

11.81.1 Constructor & Destructor Documentation

11.81.1.1 CPgaDeviceNet() CPgaDeviceNet ()

```
11.81.1.2 ~CPgaDeviceNet() ~CPgaDeviceNet ()
11.81.2 Member Function Documentation
11.81.2.1 ApplyGains() uint32_t ApplyGains ()
11.81.2.2 DefineAmplification() uint32_t DefineAmplification (
              int index,
              int amplification,
              int poti1,
              int poti2 )
11.81.2.3 DefineFrequencyRange() uint32_t DefineFrequencyRange (
              int index,
              int low,
              int high,
              int channels,
              int gain )
11.81.2.4 DefineNumAmplifications() uint32_t DefineNumAmplifications (
              int number )
\textbf{11.81.2.5} \quad \textbf{DefineNumFrequencyRanges()} \quad \texttt{uint32\_t DefineNumFrequencyRanges} \quad \textbf{(}
              int rnum )
11.81.2.6 GetAmplification() uint32_t GetAmplification (
              int index,
              [System::Runtime::InteropServices::Out] int% amplification,
              [System::Runtime::InteropServices::Out] int% poti1,
```


[System::Runtime::InteropServices::Out] int% poti2)

```
11.81.2.7 GetFrequencyRange() uint32_t GetFrequencyRange (
             int rangeIndex,
             [System::Runtime::InteropServices::Out] int% low,
             [System::Runtime::InteropServices::Out] int% high,
             [System::Runtime::InteropServices::Out] int% channels,
             [System::Runtime::InteropServices::Out] int% gain )
11.81.2.8 GetGain() uint32_t GetGain (
             int channel,
             [System::Runtime::InteropServices::Out] int% Gain,
             [System::Runtime::InteropServices::Out] int% poti1,
             [System::Runtime::InteropServices::Out] int% poti2 )
11.81.2.9 GetNumAmplifications() uint32_t GetNumAmplifications (
             [System::Runtime::InteropServices::Out] int% number )
11.81.2.10 GetNumFrequencyRanges() uint32_t GetNumFrequencyRanges (
             [System::Runtime::InteropServices::Out] int% numRanges )
11.81.2.11 SetGain() uint32_t SetGain (
             int channel,
             int Gain,
             int poti1,
             int poti2 )
```

11.82 CPositionIIDeviceNet Class Reference

CPositionIIDeviceNet is the class to control PositionII devices

Inheritance diagram for CPositionIIDeviceNet:

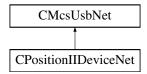

Public Member Functions

CPositionIIDeviceNet ()

Initializes a new instance of the CPositionIIDeviceNet class.

- virtual ∼CPositionIIDeviceNet ()
- !CPositionIIDeviceNet ()
- uint32_t GetCoilCommunication (uint16_t coil)

get if the communication to the coil is working

uint32 t GetOnOff (uint16 t coil)

get if the coil is switched on/off

void SwitchOnOff (uint16_t coil, uint32_t on)

switched the coild on of

uint32_t GetImplantState (uint16_t coil)

gets the implantat state

uint32_t GetImplantCurrentSetpoint (uint16_t coil)

sets the implant current setpoint

void SetImplantCurrentSetpoint (uint16 t coil, uint32 t current)

gets the implant current setpoint

uint32 t GetPowerStrength (uint16 t coil)

sets the power for the trigger pulses

void SetPowerStrength (uint16_t coil, uint32_t power)

gets the power for the trigger pulses

uint32 t GetImplantResult (uint16 t coil)

gets the last result of the implant pulse trigger

void GetRTC ([System::Runtime::InteropServices::Out]uint8_t% year, [System::Runtime::InteropServices
::Out]uint8_t% month, [System::Runtime::InteropServices::Out]uint8_t% day, [System::Runtime::Interop
Services::Out]uint8_t% hour, [System::Runtime::InteropServices::Out]uint8_t% minute, [System::Runtime
::InteropServices::Out]uint8_t% second)

Get the RTC

• void SetRTC (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)

Set the RTC

uint32_t GetStateDebugData (uint16_t coil)

get the debug queue state

void SetStateDebugData (uint16_t coil, uint32_t state)

clears/starts/stops the debug queue for a certain coil

void GetDebugData (uint16_t coil, [System::Runtime::InteropServices::Out]uint16_t% index, [System::
Runtime::InteropServices::Out]uint16_t% voltage, [System::Runtime::InteropServices::Out]uint16_t% numberofpulses, [System::Runtime::InteropServices::Out]uint16_t% mediantime)

get the oldest debug entry for a certain coil

uint32_t GetStateEventData ()

get the event queue state

void SetStateEventData (uint32_t state)

clears/starts/stops the event queue for a certain coil

get the oldest event entry

Properties

• CRFFunctionNet^ RFFunction [get]

Additional Inherited Members

11.82.1 Detailed Description

CPositionIIDeviceNet is the class to control PositionII devices

11.82.2 Constructor & Destructor Documentation

```
11.82.2.1 CPositionIIDeviceNet() CPositionIIDeviceNet ()
```

Initializes a new instance of the CPositionIIDeviceNet class.

```
11.82.2.2 ~CPositionIIDeviceNet() virtual ~CPositionIIDeviceNet () [virtual]
```

```
11.82.2.3 "!CPositionIIDeviceNet() !CPositionIIDeviceNet ( )
```

11.82.3 Member Function Documentation

```
11.82.3.1 GetCoilCommunication() uint32_t GetCoilCommunication ( uint16_t coil)
```

get if the communication to the coil is working

Parameters

coil the coil

Returns

is communicating

11.82.3.2 GetDebugData() void GetDebugData (

```
uint16_t coil,
[System::Runtime::InteropServices::Out] uint16_t% index,
[System::Runtime::InteropServices::Out] uint16_t% voltage,
[System::Runtime::InteropServices::Out] uint16_t% numberofpulses,
[System::Runtime::InteropServices::Out] uint16_t% mediantime)
```

get the oldest debug entry for a certain coil

Parameters

coil	the coil
index	the debug entry index number
voltage	the voltage applied
numberofpulses	the number of pulses detected
mediantime	the median time between pulses

11.82.3.3 GetEventData() void GetEventData (

```
[System::Runtime::InteropServices::Out] uint16_t% index,
[System::Runtime::InteropServices::Out] uint8_t% year,
[System::Runtime::InteropServices::Out] uint8_t% month,
[System::Runtime::InteropServices::Out] uint8_t% day,
[System::Runtime::InteropServices::Out] uint8_t% hour,
[System::Runtime::InteropServices::Out] uint8_t% minute,
[System::Runtime::InteropServices::Out] uint8_t% second,
[System::Runtime::InteropServices::Out] uint16_t% coil,
[System::Runtime::InteropServices::Out] uint16_t% type,
[System::Runtime::InteropServices::Out] uint16_t% value)
```

get the oldest event entry

Parameters

index	the event index number
year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second
coil	the coil
type	the event type
value	the even value

11.82.3.4 GetImplantCurrentSetpoint() uint32_t GetImplantCurrentSetpoint (uint16_t coil)

sets the implant current setpoint

Da	ro	m	Δi	ŀΔ	rc
Εи	17		ы	ш	15

coil t	he coil
--------	---------

the current

11.82.3.5 GetImplantResult() uint32_t GetImplantResult (uint16_t coil)

gets the last result of the implant pulse trigger

Parameters

```
coil the coil
```

Returns

the result

11.82.3.6 GetImplantState() uint32_t GetImplantState (uint16_t coil)

gets the implantat state

Parameters

```
coil the coil
```

Returns

get if the coil is switched on/off

Parameters

coil the coil

0 = off, 1 = on

```
11.82.3.8 GetPowerStrength() uint32_t GetPowerStrength ( uint16_t coil )
```

sets the power for the trigger pulses

Parameters

coil	the coil
COII	LITE COIL

Returns

the power

11.82.3.9 GetRTC() void GetRTC()

```
[System::Runtime::InteropServices::Out] uint8_t% year,
[System::Runtime::InteropServices::Out] uint8_t% month,
[System::Runtime::InteropServices::Out] uint8_t% day,
[System::Runtime::InteropServices::Out] uint8_t% hour,
[System::Runtime::InteropServices::Out] uint8_t% minute,
[System::Runtime::InteropServices::Out] uint8_t% second )
```

Get the RTC

Parameters

year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second

```
11.82.3.10 GetStateDebugData() uint32_t GetStateDebugData ( uint16_t coil )
```

get the debug queue state

Parameters

coil the coil

the state

```
11.82.3.11 GetStateEventData() uint32_t GetStateEventData ( )
```

get the event queue state

Returns

the state

gets the implant current setpoint

Parameters

coil	the coil
current	the current

gets the power for the trigger pulses

Parameters

coil	the coil
power	the power

Set the RTC

year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second

clears/starts/stops the debug queue for a certain coil

Parameters

coil	the coil
state	clear/start/stop

11.82.3.16 SetStateEventData() void SetStateEventData (uint32_t state)

clears/starts/stops the event queue for a certain coil

Parameters

state	clear/start/stop
-------	------------------

switched the coild on of

Parameters

coil	the coil
on	0 = off, 1 = on


11.82.4 Property Documentation

11.82.4.1 RFFunction CRFFunctionNet^ RFFunction [get]

11.83 CPositionImpDeviceNet Class Reference

CPositionImpDeviceNet is the class to access the Position/Imp devices

Inheritance diagram for CPositionImpDeviceNet:

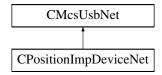

Public Member Functions

• CPositionImpDeviceNet ()

Initializes a new instance of the CPositionImpDeviceNet class.

- virtual \sim CPositionImpDeviceNet ()
- !CPositionImpDeviceNet ()
- void ConnectImp (uint32_t id)

Connect to a Imp device with a certain ID

• uint32_t ConnectedImp ()

The ID of the connected Imp device

• int32_t GetRFFrequency ()

Gets currently used RF frequency

void SetRFFrequency (int32_t frequency)

Sets the current RF frequency

• uint32_t GetDeviceList (int32_t index)

Gets the device list

• void SetDeviceList (int32_t index, uint32_t id)

Sets the device list

• uint32_t GetImpId ()

Gets the ID of the impedance measure device

void SetImpId (uint32_t id)

Sets the ID of the impedance measure device

Additional Inherited Members

11.83.1 Detailed Description

CPositionImpDeviceNet is the class to access the Position/Imp devices

11.83.2 Constructor & Destructor Documentation

```
11.83.2.1 CPositionImpDeviceNet() CPositionImpDeviceNet ()
```

Initializes a new instance of the CPositionImpDeviceNet class.

```
11.83.2.2 ~CPositionImpDeviceNet() virtual ~CPositionImpDeviceNet () [virtual]
```

```
11.83.2.3 "!CPositionImpDeviceNet() !CPositionImpDeviceNet ()
```

11.83.3 Member Function Documentation

```
11.83.3.1 ConnectedImp() uint32_t ConnectedImp ()
```

The ID of the connected Imp device

Returns

The ID

```
11.83.3.2 ConnectImp() void ConnectImp ( uint32_t id )
```

Connect to a Imp device with a certain ID

Parameters

id The ID

Gets the device list

index the index

Returns

the ID

11.83.3.4 GetImpld() uint32_t GetImpld ()

Gets the ID of the impedance measure device

Returns

the ID

11.83.3.5 GetRFFrequency() int32_t GetRFFrequency ()

Gets currently used RF frequency

Returns

The frequency

Sets the device list

Parameters

index	the index
id	the ID

```
11.83.3.7 SetImpId() void SetImpId ( uint32_t id )
```

Sets the ID of the impedance measure device

id the ID


Sets the current RF frequency

Parameters

frequency The frequency

11.84 CPPCDeviceNet Class Reference

Inheritance diagram for CPPCDeviceNet:

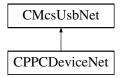

Public Member Functions

• CPPCDeviceNet (void)

Properties

- CPPCFunctionNet^ PPCFunction [get]
- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet^ McsBus_MotorControl [get]
- CMcsBus_SensorNet^ McsBus_Sensor [get]

Additional Inherited Members

11.84.1 Constructor & Destructor Documentation

11.84.1.1 CPPCDeviceNet() CPPCDeviceNet (void)

11.84.2 Property Documentation

```
11.84.2.1 McsBus CMcsBusNet^ McsBus [get]
```

11.84.2.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]


11.84.2.3 McsBus_Sensor CMcsBus_SensorNet^ McsBus_Sensor [get]

11.84.2.4 PPCFunction CPPCFunctionNet^ PPCFunction [get]

11.85 CPPCFunctionNet Class Reference

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

Inheritance diagram for CPPCFunctionNet:

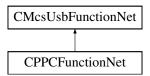

Public Member Functions

CPPCFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pPPCFunctionPointer←
 Container)

Initializes a new instance of the CPPCFunctionNet class.

- CPPCFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CPPCFunctionNet ()
- !CPPCFunctionNet ()
- int GetPumpSpeedUnit (uint16_t channel)

Reads the Pump Speed Unit

void SetPumpSpeedUnit (uint16_t channel, int SpeedUnit)

Writes the Pump Speed Unit

PP_Pump_Mode_Type_EnumNet GetPumpModeType (uint16_t channel)

Reads the Pump Mode Type.

• void SetPumpModeType (uint16_t channel, PP_Pump_Mode_Type_EnumNet PumpMode)

Writes the config string from the device.

void GetAnalogVoltageRange (uint16_t channel, [System::Runtime::InteropServices::Out]uint16_t% min_
 voltage, [System::Runtime::InteropServices::Out]uint16_t% max_voltage)

Reads the Analog Input Voltage Range

void SetAnalogVoltageRange (uint16_t channel, uint16_t min_voltage, uint16_t max_voltage)

Writes the Analog Input Voltage Range

void GetPressureRange (uint16_t channel, [System::Runtime::InteropServices::Out]int32_t% lower_
 pressure, [System::Runtime::InteropServices::Out]int32_t% upper_pressure)

Get the pressure range that is used between the analog voltage or the digital states

void SetPressureRange (uint16_t channel, int32_t lower_pressure, int32_t upper_pressure)

Get the pressure range that is used between the analog voltage or the digital states

uint16_t GetSupplyVoltage ()

Reads the current supply voltage in mV

uint16 t GetAnalogVoltage (uint16 t channel)

Reads the current analog voltage

uint16_t GetDigitalIn (uint16_t channel)

Reads the digital input state

• int GetValveActive (uint16 t valve)

Gets the valve active/inactive state

void SetValveActive (uint16_t valve, int valveActive)

Sets the valve active/inactive state

void SetPressureOffset ()

Sets the pressure offset

void LoadPressure (int32_t pressure, uint32_t options)

Loads the reservoir with a pressure

void IsBusy ([System::Runtime::InteropServices::Out]int16_t% task, [System::Runtime::InteropServices::
 Out]int16_t% wait)

Is the PPC busy with a task

void FirePressurePulse (int32_t duration, int32_t nextpressure)

Fire a pressure pulse from the reservoir

int32_t MeasureReservoir ()

Measures the reservoir pressure

Additional Inherited Members

11.85.1 Detailed Description

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

11.85.2 Constructor & Destructor Documentation

Initializes a new instance of the CPPCFunctionNet class.

11.85.3 Member Function Documentation

Fire a pressure pulse from the reservoir

Parameters

duration	The pulse duration (valves open)
nextpressure	The next pressure

Reads the current analog voltage

Parameters

Returns

The Analog Voltage

[System::Runtime::InteropServices::Out] uint16_t% min_voltage, [System::Runtime::InteropServices::Out] uint16_t% max_voltage)

Reads the Analog Input Voltage Range

Parameters

channel	The Channel Number
min_voltage	The voltage that should be seen as the minimum voltage
max_voltage	The voltage that should be seen as the maximum voltage

Reads the digital input state

Parameters

channel The Channel Number

Returns

The Digital State

Get the pressure range that is used between the analog voltage or the digital states

Parameters

channel	The Channel Number
lower_pressure	The lower border of the pressure range
upper_pressure	The upper border of the pressure range

```
11.85.3.6 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType ( uint16_t channel )
```

Reads the Pump Mode Type.

Parameters

channel	The Channel Number

Returns

The Pump Mode Type.

Reads the Pump Speed Unit

Parameters

Returns

The Speed Unit

11.85.3.8 GetSupplyVoltage() uint16_t GetSupplyVoltage ()

Reads the current supply voltage in mV

Returns

The supply voltage

```
11.85.3.9 GetValveActive() int GetValveActive ( uint16_t valve )
```

Gets the valve active/inactive state

Parameters

```
valve The valve number
```

Returns

The valve state

Is the PPC busy with a task

Parameters

task	The task state
wait	The wait state

Loads the reservoir with a pressure

Parameters

pressure	The pressure
options	The options: end with 0=regulate on patch 1=regulate on reservoir

11.85.3.12 MeasureReservoir() int32_t MeasureReservoir ()

Measures the reservoir pressure

Returns

The pressure

Writes the Analog Input Voltage Range

Parameters

channel	The Channel Number
min_voltage	The voltage that should be seen as the minimum voltage
max_voltage	The voltage that should be seen as the maximum voltage

11.85.3.14 SetPressureOffset() void SetPressureOffset ()

Sets the pressure offset

Get the pressure range that is used between the analog voltage or the digital states

Parameters

channel	The Channel Number
lower_pressure	The lower border of the pressure range
upper_pressure	The upper border of the pressure range

Writes the config string from the device.

Parameters

channel	The Channel Number
PumpMode	The Pump Mode Type.

Writes the Pump Speed Unit

Parameters

channel	The Channel Number
SpeedUnit	The Speed Unit


Sets the valve active/inactive state

Parameters

valve	The valve number
valveActive	The valve state

11.86 CPPS_DeviceNet Class Reference

Inheritance diagram for CPPS_DeviceNet:

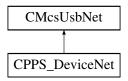

Public Member Functions

CPPS_DeviceNet (void)

Properties

- CPPS_FunctionNet^ PPS_Function [get]
- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet^ McsBus_MotorControl [get]
- CMcsBus_SensorNet^ McsBus_Sensor [get]

Additional Inherited Members

11.86.1 Constructor & Destructor Documentation

```
11.86.1.1 CPPS_DeviceNet() CPPS_DeviceNet (
void )
```

11.86.2 Property Documentation

11.86.2.1 McsBus CMcsBusNet^ McsBus [get]


11.86.2.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.86.2.3 McsBus_Sensor CMcsBus_SensorNet^ McsBus_Sensor [get]

11.86.2.4 PPS_Function CPPS_FunctionNet^ PPS_Function [get]

11.87 CPPS_FunctionNet Class Reference

Inheritance diagram for CPPS_FunctionNet:

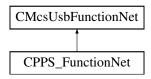

Public Member Functions

- CPPS_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] cPPS_FunctionPointer←
 Container)
- CPPS FunctionNet (CMcsUsbNet[^] mcsusb)
- · void SetPumpMaxSpeed (unsigned short index, unsigned short maxspeed)
- unsigned short GetPumpMaxSpeed (unsigned short index)
- void SetPumpSpeedUnit (unsigned short index, int speedunit)
- int GetPumpSpeedUnit (unsigned short index)
- void SetPumpModeType (unsigned short index, PP_Pump_Mode_Type_EnumNet type)
- PP_Pump_Mode_Type_EnumNet GetPumpModeType (unsigned short index)
- void SetPumpCouple (unsigned int i)
- unsigned int GetPumpCouple ()
- void SetPumpEnableSpeedRatio (unsigned int enable)
- unsigned int GetPumpEnableSpeedRatio ()
- · void SetPumpManualOnOff (unsigned short index, unsigned int onoff)
- unsigned int GetPumpManualOnOff (unsigned short index)
- void SetPumpFunctionSpeeds (unsigned short index, short offspeed, short onspeed)
- void GetPumpFunctionSpeeds (unsigned short index, [System::Runtime::InteropServices::Out]short% offspeed, [System::Runtime::InteropServices::Out]short% onspeed)
- void SetPumpSpeedRatio (int ratio)
- int GetPumpSpeedRatio ()
- void SetPumpFastOnOff (unsigned short index, unsigned int onoff)
- unsigned int GetPumpFastOnOff (unsigned short index)
- void SetPumpFastSpeed (unsigned short index, short fastspeed)
- short GetPumpFastSpeed (unsigned short index)
- void SetAnalogVoltages (unsigned short index, unsigned short minvoltage, unsigned short maxvoltage)
- void GetAnalogVoltages (unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% minvoltage, [System::Runtime::InteropServices::Out]unsigned short% maxvoltage)
- void SetUseBubble (unsigned short index, unsigned int usebubble)
- unsigned int GetUseBubble (unsigned short index)
- unsigned short GetSupplyVoltage ()
- unsigned short GetAnalogVoltage (unsigned short index)
- unsigned short GetDigitalIn (unsigned short index)
- unsigned short GetBubbleState ()

Additional Inherited Members

11.87.1 Constructor & Destructor Documentation

```
11.87.1.1 CPPS_FunctionNet() [1/2] CPPS_FunctionNet (
             CMcsUsbNet^ mcsusb,
             {\tt CMcsUsbFunctionPointerContainer}^{\wedge} \ \ cPPS\_FunctionPointerContainer \ )
11.87.1.2 CPPS_FunctionNet() [2/2] CPPS_FunctionNet (
             CMcsUsbNet^ mcsusb )
11.87.2 Member Function Documentation
11.87.2.1 GetAnalogVoltage() unsigned short GetAnalogVoltage (
             unsigned short index)
11.87.2.2 GetAnalogVoltages() void GetAnalogVoltages (
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned short% minvoltage,
             [System::Runtime::InteropServices::Out] unsigned short% maxvoltage )
11.87.2.3 GetBubbleState() unsigned short GetBubbleState ( )
11.87.2.4 GetDigitalIn() unsigned short GetDigitalIn (
             unsigned short index )
11.87.2.5 GetPumpCouple() unsigned int GetPumpCouple ()
11.87.2.6 GetPumpEnableSpeedRatio() unsigned int GetPumpEnableSpeedRatio ( )
11.87.2.7 GetPumpFastOnOff() unsigned int GetPumpFastOnOff (
             unsigned short index)
```


```
11.87.2.8 GetPumpFastSpeed() short GetPumpFastSpeed (
             unsigned short index)
11.87.2.9 GetPumpFunctionSpeeds() void GetPumpFunctionSpeeds (
             unsigned short index,
             [System::Runtime::InteropServices::Out] short% offspeed,
             [System::Runtime::InteropServices::Out] short% onspeed )
\textbf{11.87.2.10} \quad \textbf{GetPumpManualOnOff()} \quad \texttt{unsigned int GetPumpManualOnOff} \quad \textbf{(}
             unsigned short index)
11.87.2.11 GetPumpMaxSpeed() unsigned short GetPumpMaxSpeed (
             unsigned short index)
11.87.2.12 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType (
             unsigned short index)
11.87.2.13 GetPumpSpeedRatio() int GetPumpSpeedRatio ( )
11.87.2.14 GetPumpSpeedUnit() int GetPumpSpeedUnit (
             unsigned short index)
11.87.2.15 GetSupplyVoltage() unsigned short GetSupplyVoltage ( )
11.87.2.16 GetUseBubble() unsigned int GetUseBubble (
             unsigned short index)
```

```
\textbf{11.87.2.17} \quad \textbf{SetAnalogVoltages()} \quad \texttt{void SetAnalogVoltages ()}
              unsigned short index,
              unsigned short minvoltage,
              unsigned short maxvoltage )
11.87.2.18 SetPumpCouple() void SetPumpCouple (
              unsigned int i )
11.87.2.19 SetPumpEnableSpeedRatio() void SetPumpEnableSpeedRatio (
              unsigned int enable )
11.87.2.20 SetPumpFastOnOff() void SetPumpFastOnOff (
              unsigned short index,
              unsigned int onoff )
\textbf{11.87.2.21} \quad \textbf{SetPumpFastSpeed()} \quad \texttt{void} \ \texttt{SetPumpFastSpeed} \ \ \textbf{(}
              unsigned short index,
              short fastspeed )
11.87.2.22 SetPumpFunctionSpeeds() void SetPumpFunctionSpeeds (
              unsigned short index,
              short offspeed,
              short onspeed )
11.87.2.23 SetPumpManualOnOff() void SetPumpManualOnOff (
              unsigned short index,
              unsigned int onoff )
11.87.2.24 SetPumpMaxSpeed() void SetPumpMaxSpeed (
              unsigned short index,
              unsigned short {\it maxspeed} )
```

11.88 CPPSDeviceNet Class Reference

CPPS4plus1DeviceNet is the to control the MCS HLA device

Inheritance diagram for CPPSDeviceNet:

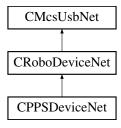

Public Member Functions

• CPPSDeviceNet (void)

Additional Inherited Members

11.88.1 Detailed Description

CPPS4plus1DeviceNet is the to control the MCS HLA device


11.88.2 Constructor & Destructor Documentation

```
11.88.2.1 CPPSDeviceNet() CPPSDeviceNet (
void )
```

11.89 CProgramPressureCurveNet Class Reference

CProgramPressureCurveNet is the class to program pressure curves

Inheritance diagram for CProgramPressureCurveNet:

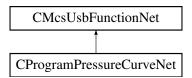

Public Member Functions

- CProgramPressureCurveNet (CMcsUsbNet[^] mcsusb)
 Initializes a new instance of the CPulseGeneratorFunctionNet class.
- virtual ~CProgramPressureCurveNet (void)
- !CProgramPressureCurveNet (void)
- void Program (unsigned char busnumber, unsigned char busaddress, int32_t channel, array< int32_t $>^{\land}$ pressures, array< int32_t $>^{\land}$ durations)
- void SetRepeats (unsigned char busnumber, unsigned char busaddress, int32_t channel, uint32_t repeats)
- unsigned int GetRepeats (unsigned char busnumber, unsigned char busaddress, int32 t channel)

Additional Inherited Members

11.89.1 Detailed Description

CProgramPressureCurveNet is the class to program pressure curves

11.89.2 Constructor & Destructor Documentation

```
11.89.2.1 CProgramPressureCurveNet() CProgramPressureCurveNet (
CMcsUsbNet^ mcsusb )
```


Initializes a new instance of the CPulseGeneratorFunctionNet class.

```
11.89.2.2 ~CProgramPressureCurveNet() virtual ~CProgramPressureCurveNet (
             void ) [virtual]
11.89.2.3 "!CProgramPressureCurveNet() !CProgramPressureCurveNet (
             void )
11.89.3 Member Function Documentation
11.89.3.1 GetRepeats() unsigned int GetRepeats (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel )
11.89.3.2 Program() void Program (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel,
             array< int32_t >^{\wedge} pressures,
             array< int32_t >^{\wedge} steps,
             array< int16_t >^{\land} durations )
11.89.3.3 SetRepeats() void SetRepeats (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel,
             uint32_t repeats )
```

11.90 CPulseGeneratorFunctionNet Class Reference

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

Inheritance diagram for CPulseGeneratorFunctionNet:

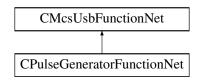

Public Member Functions

 CPulseGeneratorFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pPulse← GeneratorFunctionPointerContainer)

Initializes a new instance of the CPulseGeneratorFunctionNet class.

- CPulseGeneratorFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CPulseGeneratorFunctionNet ()
- !CPulseGeneratorFunctionNet ()
- int32_t GetPeriod (int32_t generator_number)

Reads the generator period

• void SetPeriod (int32_t generator_number, int32_t period_in_samples)

Writes the generator period

int32_t GetPulseLength (int32_t generator_number)

Reads the generator pulse length

void SetPulseLength (int32 t generator number, int32 t pulselength in 10us)

Writes the generator pulse length

void GetModeSelect (int32_t generator_number, [System::Runtime::InteropServices::Out]PulseGenerator
 —Mode_EnumNet% mode, [System::Runtime::InteropServices::Out]int32_t% digitalchannel)

Reads the generator mode

void SetModeSelect (int32_t generator_number, PulseGenerator_Mode_EnumNet mode, int32_t digitalchannel)

Writes the generator mode

Additional Inherited Members

11.90.1 Detailed Description

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

11.90.2 Constructor & Destructor Documentation

```
11.90.2.1 CPulseGeneratorFunctionNet() [1/2] CPulseGeneratorFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pPulseGeneratorFunctionPointerContainer)
```

Initializes a new instance of the CPulseGeneratorFunctionNet class.

```
11.90.2.2 CPulseGeneratorFunctionNet() [2/2] CPulseGeneratorFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.90.2.3 ~CPulseGeneratorFunctionNet() virtual ~CPulseGeneratorFunctionNet ( ) [virtual]
```

11.90.2.4 "!CPulseGeneratorFunctionNet() !CPulseGeneratorFunctionNet ()

11.90.3 Member Function Documentation

Reads the generator mode

Parameters

generator_number	The generator number
mode	The generator mode
digitalchannel	The digital in channel used as gate

Reads the generator period

Parameters

generator_number	The generator number
------------------	----------------------

Returns

The period

Reads the generator pulse length

Parameters

Returns

The pulse length

```
11.90.3.4 SetModeSelect() void SetModeSelect (
    int32_t generator_number,
    PulseGenerator_Mode_EnumNet mode,
    int32_t digitalchannel )
```

Writes the generator mode

Parameters

generator_number	The generator number
mode	The generator mode
digitalchannel	The digital in channel used as gate

Writes the generator period

Parameters

generator_number	The generator number
period_in_samples	The period


Writes the generator pulse length

Parameters

generator_number	The generator number
pulselength_in_10us	The pulse length

11.91 CRadioControledDevicesNet Class Reference

Inheritance diagram for CRadioControledDevicesNet:

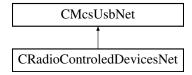

Public Member Functions

- CRadioControledDevicesNet (void)
- bool HasRadioControl ()
- array< unsigned short > ^ GetDeviceNames ()
- void ConnectDevice (unsigned short sn)
- void DisConnectDevice ()
- bool StillConnected ()
- void SetFrequency (unsigned short frequency)
- unsigned short GetFrequency ()

Protected Member Functions

• CRadioControledDevicesNet (CRadioControledDevices *pRadioControled)

Additional Inherited Members

11.91.1 Constructor & Destructor Documentation

```
11.91.1.1 CRadioControledDevicesNet() [1/2] CRadioControledDevicesNet ( void )
```

```
11.91.1.2 CRadioControledDevicesNet() [2/2] CRadioControledDevicesNet (
CRadioControledDevices * pRadioControled) [protected]
```

11.91.2 Member Function Documentation

```
11.91.2.1 ConnectDevice() void ConnectDevice (
unsigned short sn )
```

11.91.2.2 DisConnectDevice() void DisConnectDevice ()

11.92 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference

Public Member Functions

- CRegionOfInterestRect (int left, int top, int right, int bottom)
- CRegionOfInterestRect ^ DeepCopy ()

Public Attributes

```
• int m_Left
```

- int m_Top
- int m_Right
- int m_Bottom

11.92.1 Constructor & Destructor Documentation

```
11.92.1.1 CRegionOfInterestRect() CRegionOfInterestRect (
          int left,
          int top,
          int right,
          int bottom )
```

11.92.2 Member Function Documentation

11.92.2.1 DeepCopy() CRegionOfInterestRect ^ DeepCopy ()

11.92.3 Member Data Documentation

11.92.3.1 m_Bottom int m_Bottom


11.92.3.2 m_Left int m_Left

11.92.3.3 m_Right int m_Right

11.92.3.4 m_Top int m_Top

11.93 CRetinaLedDeviceNet Class Reference

Inheritance diagram for CRetinaLedDeviceNet:

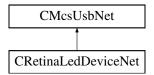

Public Member Functions

- CRetinaLedDeviceNet ()
- ∼CRetinaLedDeviceNet ()
- unsigned int SetTrigger (int enable)
- unsigned int SetLED (unsigned long long pattern)
- unsigned int SetTablepointer (int position)
- unsigned int GetTablepointer (int % position)
- unsigned int ClearTable ()
- unsigned int AddTableEntry (unsigned long long pattern)
- unsigned int AddLoopEntry (unsigned short repeats, unsigned short steps_back)
- unsigned int SetRepeat (int repeat)
- unsigned int SetLumi (int lumi)
- unsigned int SetPersistency (unsigned int persistency)

Additional Inherited Members


```
11.93.1 Constructor & Destructor Documentation
```

```
11.93.1.1 CRetinaLedDeviceNet() CRetinaLedDeviceNet ()
11.93.1.2 \sim CRetinaLedDeviceNet() \sim CRetinaLedDeviceNet ()
11.93.2 Member Function Documentation
11.93.2.1 AddLoopEntry() unsigned int AddLoopEntry (
             unsigned short repeats,
             unsigned short steps_back )
11.93.2.2 AddTableEntry() unsigned int AddTableEntry (
             unsigned long long pattern )
11.93.2.3 ClearTable() unsigned int ClearTable ( )
11.93.2.4 GetTablepointer() unsigned int GetTablepointer (
             int % position )
11.93.2.5 SetLED() unsigned int SetLED (
             unsigned long long pattern )
11.93.2.6 SetLumi() unsigned int SetLumi (
             int lumi )
```

11.94 CRFFunctionNet Class Reference

CRFFunctionNet is the class to control RF devices

Inheritance diagram for CRFFunctionNet:

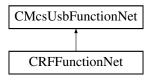

Public Member Functions

CRFFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pRFFunctionPointer←
 Container)

Initializes a new instance of the CRFFunctionNet class.

- CRFFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CRFFunctionNet ()
- !CRFFunctionNet ()
- uint32_t GetBaseFrequency (CFirmwareDestinationNet destination)

gets the base advertise frequency

void SetBaseFrequency (CFirmwareDestinationNet destination, uint32_t frequency)

gets the base advertise frequency

• uint32 t GetWorkingFrequency ()

gets the working frequency

void SetWorkingFrequency (uint32_t frequency)

sets the working frequency

array< uint32_t > ^ GetAvailableDeviceListEx (int list_Length)

get a list of available devices

array< uint32_t > ^ GetAvailableDeviceList ()

```
get a list of available devices
    - array< uint32_t > ^{\wedge} GetAvailableStateListEx (int list_Length)
          get a list of the states of the available devices

    array< uint32_t > ^ GetAvailableStateList ()

          get a list of the states of the available devices
    • void Connect (uint32_t sn)
          connect to a RF device, use 0 to disconnect

    uint32_t GetConnectedDevice ()

          get connect RF device, 0 = no device connected
    • uint32_t GetState ()
          get connection state

    void SetTestMode (uint32_t mode)

          set test mode

    uint32_t GetTestMode ()

          gets test mode
Additional Inherited Members
11.94.1 Detailed Description
CRFFunctionNet is the class to control RF devices
11.94.2 Constructor & Destructor Documentation
```

```
11.94.2.1 CRFFunctionNet() [1/2] CRFFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pRFFunctionPointerContainer)
```

Initializes a new instance of the CRFFunctionNet class.

11.94.3 Member Function Documentation

```
11.94.3.1 Connect() void Connect ( uint32\_t sn )
```

connect to a RF device, use 0 to disconnect

Parameters sn
11.94.3.2 GetAvailableDeviceList() array <uint32_t> ^ GetAvailableDeviceList ()</uint32_t>
get a list of available devices
Returns
<pre>11.94.3.3 GetAvailableDeviceListEx() array<uint32_t> ^ GetAvailableDeviceListEx (int list_Length)</uint32_t></pre>
get a list of available devices
Parameters
list_Length The maximal length of list.
Returns
11.94.3.4 GetAvailableStateList() array <uint32_t> ^ GetAvailableStateList ()</uint32_t>
get a list of the states of the available devices
Returns
11.94.3.5 GetAvailableStateListEx() array <uint32_t> ^ GetAvailableStateListEx (</uint32_t>

Generated by Doxygen

int list_Length)

get a list of the states of the available devices

Pa	ra	m	۵	łα	re
Гα	ıα		е	ш	15

list_Length	The maximal length of list.	
-------------	-----------------------------	--

Returns

```
11.94.3.6 GetBaseFrequency() uint32_t GetBaseFrequency ( CFirmwareDestinationNet destination)
```

gets the base advertise frequency

Parameters

destination

Returns

11.94.3.7 GetConnectedDevice() uint32_t GetConnectedDevice ()

get connect RF device, 0 = no device connected

Returns

11.94.3.8 GetState() uint32_t GetState ()

get connection state

Returns

```
11.94.3.9 GetTestMode() uint32_t GetTestMode ( )
gets test mode
Returns
11.94.3.10 GetWorkingFrequency() uint32_t GetWorkingFrequency ( )
gets the working frequency
Returns
11.94.3.11 SetBaseFrequency() void SetBaseFrequency (
             CFirmwareDestinationNet destination,
             uint32_t frequency )
gets the base advertise frequency
Parameters
 destination
 frequency
11.94.3.12 SetTestMode() void SetTestMode (
             uint32_t mode )
set test mode
Parameters
 mode
11.94.3.13 SetWorkingFrequency() void SetWorkingFrequency (
```

sets the working frequency


uint32_t frequency)

Parameters

frequency

11.95 CRobo_FYIProgram_FunctionNet Class Reference

Inheritance diagram for CRobo_FYIProgram_FunctionNet:

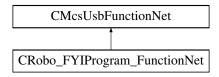

Public Member Functions

- CRobo_FYIProgram_FunctionNet (CMcsUsbNet[^] mcsusb)
- void SetValve1 (unsigned char index, unsigned int valve1)
- unsigned int GetValve1 (unsigned char index)
- void SetValve2 (unsigned char index, unsigned int valve2)
- unsigned int GetValve2 (unsigned char index)
- · void SetLength (unsigned char index, int length)
- int GetLength (unsigned char index)
- void Start ()
- · int GetState ()

Additional Inherited Members

11.95.1 Constructor & Destructor Documentation


```
11.95.1.2 CRobo_FYIProgram_FunctionNet() [2/2] CRobo_FYIProgram_FunctionNet (
CMcsUsbNet^ mcsusb )
```

11.95.2 Member Function Documentation

```
\textbf{11.95.2.1} \quad \textbf{GetLength()} \quad \texttt{int GetLength ()}
              unsigned char index )
11.95.2.2 GetState() int GetState ()
11.95.2.3 GetValve1() unsigned int GetValve1 (
              unsigned char index)
11.95.2.4 GetValve2() unsigned int GetValve2 (
              unsigned char index)
11.95.2.5 SetLength() void SetLength (
              unsigned char index,
              int length )
11.95.2.6 SetValve1() void SetValve1 (
              unsigned char index,
              unsigned int valve1 )
11.95.2.7 SetValve2() void SetValve2 (
              unsigned char index,
              unsigned int valve2)
11.95.2.8 Start() void Start ()
```

11.96 CRobo_FYITemp_FunctionNet Class Reference

Inheritance diagram for CRobo_FYITemp_FunctionNet:

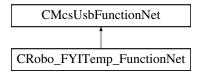

Public Member Functions

- CRobo_FYITemp_FunctionNet (CMcsUsbNet[^] mcsusb)
- void SetRegulatorOnOff (unsigned char index, int onoff)
- int GetRegulatorOnOff (unsigned char index)
- void SetSollTemp (unsigned char index, int temp)
- int GetSollTemp (unsigned char index)
- void SetPCoeff (unsigned char index, int pcoeff)
- int GetPCoeff (unsigned char index)
- void SetICoeff (unsigned char index, int icoeff)
- int GetlCoeff (unsigned char index)
- void SetMaxPower (unsigned char index, int power)
- int GetMaxPower (unsigned char index)

Additional Inherited Members

11.96.1 Constructor & Destructor Documentation

```
11.96.1.1 CRobo_FYITemp_FunctionNet() CRobo_FYITemp_FunctionNet (
CMcsUsbNet^ mcsusb)
```

11.96.2 Member Function Documentation

```
11.96.2.1 GetlCoeff() int GetICoeff (
          unsigned char index )
```

```
11.96.2.2 GetMaxPower() int GetMaxPower ( unsigned char index )
```


```
11.96.2.3 GetPCoeff() int GetPCoeff (
unsigned char index)
```

```
11.96.2.4 GetRegulatorOnOff() int GetRegulatorOnOff (
    unsigned char index )
```

```
11.96.2.5 GetSollTemp() int GetSollTemp (
             unsigned char index )
11.96.2.6 SetlCoeff() void SetICoeff (
             unsigned char index,
             int icoeff )
11.96.2.7 SetMaxPower() void SetMaxPower (
             unsigned char index,
             int power )
11.96.2.8 SetPCoeff() void SetPCoeff (
             unsigned char index,
             int pcoeff )
11.96.2.9 SetRegulatorOnOff() void SetRegulatorOnOff (
             unsigned char index,
             int onoff )
11.96.2.10 SetSollTemp() void SetSollTemp (
             unsigned char index,
             int temp )
```

11.97 CRoboDacqNet Class Reference

Inheritance diagram for CRoboDacqNet:

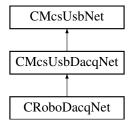

Public Member Functions

- CRoboDacqNet (void)
- CRoboDacqNet (CRoboDeviceNet[^] robodevice)
- void RunTable ()
- void RunTable (int timeout)
- void StopTable ()
- void StopTable (int timeout)
- void CancelTableLoop ()
- void CancelTableLoopAndStopTable ()
- void SetConfigurationBit (unsigned short bit, bool value)
- void SetConfigurationBitSupply (bool value)
- void SetConfigurationBitRelais (bool value)
- void SetConfigurationBitStream (bool value)
- void SetConfigurationBitAxc (bool value)
- void SetConfigurationBitCC_Gen (bool value)
- void SetConfigurationBitCV Gen (bool value)
- void SetConfigurationBitRC_Gen (bool value)
- void SetConfigurationBitRV Gen (bool value)
- void SetConfigurationBitBlu_Led (bool value)
- void SetConfigurationBitRed Led (bool value)
- void SetConfigurationBitBlu_LedToggleSlow (bool value)
- void SetConfigurationBitRed LedToggleSlow (bool value)
- void SetConfigurationBitRlu_LedToggleFast (bool value)
- void SetConfigurationBitRed_LedToggleFast (bool value)
- void SetConfigurationBitRed_LedSaturation (bool value)
- void SetSimulation (unsigned int enable)
- void SetUClamp (int uClamp)
- void SetIClamp (int iClamp)
- void SetPGain (int pGain)
- void SetlGain (int iGain)
- void SetFilter (int filter)
- void SetUVOffset (int UVOffset)
- void SetUCOffset (int UCOffset)
- void SetICOffset (int ICOffset)
- void SetCrossTalkOffset (int CrossTalk)
- void SetXGain (int xGain)
- void SetCrossTalkOptimum (int cxOptimum)
- void SetRecordingNumber (unsigned int recordingNumber)
- void ClampAmpRestart ()
- void DoRamp (int startValue, int endValue, int duration, int mode)
- unsigned int GetClampAmpSerialNumber ()
- unsigned int GetConfigurationBits ()
- bool GetConfigurationBit (unsigned short bit)
- bool GetConfigurationBitSupply ()
- bool GetConfigurationBitRelais ()
- bool GetConfigurationBitStream ()
- bool GetConfigurationBitAxc ()
- · bool GetConfigurationBitCC Gen ()
- bool GetConfigurationBitCV_Gen ()
- · bool GetConfigurationBitRC Gen ()
- bool GetConfigurationBitRV Gen ()
- bool GetConfigurationBitBlu Led ()
- bool GetConfigurationBitRed_Led ()
- bool GetConfigurationBitBlu_LedToggleSlow ()

- bool GetConfigurationBitRed_LedToggleSlow ()
- bool GetConfigurationBitBlu_LedToggleFast ()
- bool GetConfigurationBitRed LedToggleFast ()
- bool GetConfigurationBitRed_LedSaturation ()
- unsigned int GetSimulation ()
- int GetUClamp ()
- int GetlClamp ()
- int GetPGain ()
- int GetlGain ()
- int GetFilter ()
- int GetUVOffset ()
- int GetUCOffset ()
- int GetICOffset ()
- int GetCrossTalkOffset ()
- int GetXGain ()
- int GetCrossTalkOptimum ()
- unsigned int GetRecordingNumber ()
- int GetResistanceC ()
- int GetResistanceV ()
- int GetCapacityC ()
- · int GetCapacityV ()
- int GetCapacityX ()
- int GetUV ()
- int GetUC ()
- int GetIC ()
- int GetNUV_MS ()
- int GetNUC MS ()
- int GetNIC MS ()
- void SetAllDigout (uint32_t value)
- uint32_t GetAllDigout ()
- · void SetCommand (unsigned char command, int value)
- int GetCommand (unsigned char command)
- void SetDigout (uint16_t index, bool enable)
- bool GetDigout (uint16 t index)
- · void TableDefBegin ()
- void TableDefEnd ()
- void Table_Wait (unsigned int tableWait)
- void SetDownsampleFactor (int index, int downsample_factor)
- void SetFilterCoeffs (int index, array< int >^ coeffs)
- void SetNoFilterCoeffs (int index)
- int GetDownsampleFactor (int index)
- array< int > ^ GetFilterCoeffs (int index)
- void Emu_SetElectrodeResists (int emuElectrodeResist)
- void Emu SetCellResists (int emuCellResist)
- void Emu_SetCellCapacity (int emuCellCapacity)
- void Emu_SetCellPotential (int emuCellPotential)
- void Emu_SetNoise (int emuNoise)
- int Emu_GetElectrodeResists ()
- int Emu_GetCellResists ()
- int Emu_GetCellCapacity ()int Emu_GetCellPotential ()
- int Emu_GetNoise ()
- void SetDisplayText (int index, String[^] displayText)
- void SetScreen (int screen)
- void UpdateDisplay ()
- String \(^\) GetDisplayText (int index)
- int GetScreen ()
- int GetUpdateDisplay ()

Additional Inherited Members

11.97.1 Constructor & Destructor Documentation

```
11.97.1.1 CRoboDacqNet() [1/2] CRoboDacqNet (
             void )
11.97.1.2 CRoboDacqNet() [2/2] CRoboDacqNet (
             CRoboDeviceNet^ robodevice )
11.97.2 Member Function Documentation
11.97.2.1 CancelTableLoop() void CancelTableLoop ( )
11.97.2.2 CancelTableLoopAndStopTable() void CancelTableLoopAndStopTable ( )
11.97.2.3 ClampAmpRestart() void ClampAmpRestart ( )
\textbf{11.97.2.4} \quad \textbf{DoRamp()} \quad \texttt{void DoRamp (}
             int startValue,
             int endValue,
             int duration,
             int mode )
11.97.2.5 Emu_GetCellCapacity() int Emu_GetCellCapacity ( )
11.97.2.6 Emu_GetCellPotential() int Emu_GetCellPotential ( )
```

```
11.97.2.7 Emu_GetCellResists() int Emu_GetCellResists ( )
11.97.2.8 Emu_GetElectrodeResists() int Emu_GetElectrodeResists ( )
11.97.2.9 Emu_GetNoise() int Emu_GetNoise ( )
11.97.2.10 Emu_SetCellCapacity() void Emu_SetCellCapacity (
             int emuCellCapacity )
11.97.2.11 Emu_SetCellPotential() void Emu_SetCellPotential (
            int emuCellPotential )
11.97.2.12 Emu_SetCellResists() void Emu_SetCellResists (
             int emuCellResist )
11.97.2.13 Emu_SetElectrodeResists() void Emu_SetElectrodeResists (
             int emuElectrodeResist )
11.97.2.14 Emu_SetNoise() void Emu_SetNoise (
            int emuNoise )
11.97.2.15 GetAllDigout() uint32_t GetAllDigout ( )
11.97.2.16 GetCapacityC() int GetCapacityC ( )
```

```
11.97.2.17 GetCapacityV() int GetCapacityV ( )
11.97.2.18 GetCapacityX() int GetCapacityX ( )
11.97.2.19 GetClampAmpSerialNumber() unsigned int GetClampAmpSerialNumber ( )
11.97.2.20 GetCommand() int GetCommand (
                                                     unsigned char command )
11.97.2.21 GetConfigurationBit() bool GetConfigurationBit (
                                                     unsigned short bit )
\textbf{11.97.2.22} \quad \textbf{GetConfigurationBitAxc()} \quad \texttt{bool GetConfigurationBitAxc ()}
11.97.2.23 GetConfigurationBitBlu_Led() bool GetConfigurationBitBlu_Led ( )
\textbf{11.97.2.24} \quad \textbf{GetConfigurationBitBlu\_LedToggleFast()} \quad \texttt{bool GetConfigurationBitBlu\_LedToggleFast ()} \quad \texttt{bool BitBlu\_LedToggleFast ()} \quad \texttt{bool BitBlu\_
11.97.2.25 GetConfigurationBitBlu_LedToggleSlow() bool GetConfigurationBitBlu_LedToggleSlow ( )
11.97.2.26 GetConfigurationBitCC_Gen() bool GetConfigurationBitCC_Gen ( )
11.97.2.27 GetConfigurationBitCV_Gen() bool GetConfigurationBitCV_Gen ( )
```

11.97.2.28	<pre>GetConfigurationBitRC_Gen() bool GetConfigurationBitRC_Gen ()</pre>
11.97.2.29	GetConfigurationBitRed_Led() bool GetConfigurationBitRed_Led ()
11.97.2.30	GetConfigurationBitRed_LedSaturation() bool GetConfigurationBitRed_LedSaturation ()
11.97.2.31	GetConfigurationBitRed_LedToggleFast() bool GetConfigurationBitRed_LedToggleFast ()
11.97.2.32	GetConfigurationBitRed_LedToggleSlow() bool GetConfigurationBitRed_LedToggleSlow ()
11.97.2.33	GetConfigurationBitRelais() bool GetConfigurationBitRelais ()
11.97.2.34	GetConfigurationBitRV_Gen() bool GetConfigurationBitRV_Gen ()
11.97.2.35	GetConfigurationBits() unsigned int GetConfigurationBits ()
11.97.2.36	GetConfigurationBitStream() bool GetConfigurationBitStream ()
11.97.2.37	GetConfigurationBitSupply() bool GetConfigurationBitSupply ()
11.97.2.38	<pre>GetCrossTalkOffset() int GetCrossTalkOffset ()</pre>

```
11.97.2.39 GetCrossTalkOptimum() int GetCrossTalkOptimum ( )
\textbf{11.97.2.40} \quad \textbf{GetDigout()} \quad \texttt{bool GetDigout} \quad (
              uint16_t index )
11.97.2.41 GetDisplayText() String ^ GetDisplayText (
              int index )
11.97.2.42 GetDownsampleFactor() int GetDownsampleFactor (
              int index )
11.97.2.43 GetFilter() int GetFilter ( )
11.97.2.44 GetFilterCoeffs() array<int> ^{\wedge} GetFilterCoeffs (
              int index )
11.97.2.45 GetIC() int GetIC ( )
11.97.2.46 GetlClamp() int GetIClamp ()
11.97.2.47 GetICOffset() int GetICOffset ( )
11.97.2.48 GetlGain() int GetlGain ()
11.97.2.49 GetNIC_MS() int GetNIC_MS ()
```

```
\textbf{11.97.2.50} \quad \textbf{GetNUC\_MS()} \quad \texttt{int GetNUC\_MS ()}
11.97.2.51 GetNUV_MS() int GetNUV_MS ()
11.97.2.52 GetPGain() int GetPGain ( )
11.97.2.53 GetRecordingNumber() unsigned int GetRecordingNumber ()
11.97.2.54 GetResistanceC() int GetResistanceC ( )
11.97.2.55 GetResistanceV() int GetResistanceV ()
11.97.2.56 GetScreen() int GetScreen ()
11.97.2.57 GetSimulation() unsigned int GetSimulation ()
11.97.2.58 GetUC() int GetUC ( )
11.97.2.59 GetUClamp() int GetUClamp ()
11.97.2.60 GetUCOffset() int GetUCOffset ()
```

```
11.97.2.61 GetUpdateDisplay() int GetUpdateDisplay ( )
11.97.2.62 GetUV() int GetUV ( )
11.97.2.63 GetUVOffset() int GetUVOffset ( )
11.97.2.64 GetXGain() int GetXGain ( )
11.97.2.65 RunTable() [1/2] void RunTable ( )
11.97.2.66 RunTable() [2/2] void RunTable (
            int timeout )
11.97.2.67 SetAllDigout() void SetAllDigout (
            uint32_t value )
11.97.2.68 SetCommand() void SetCommand (
             unsigned char command,
             int value )
11.97.2.69 SetConfigurationBit() void SetConfigurationBit (
            unsigned short bit,
             bool value )
11.97.2.70 SetConfigurationBitAxc() void SetConfigurationBitAxc (
             bool value )
```

```
11.97.2.71 SetConfigurationBitBlu_Led() void SetConfigurationBitBlu_Led (
                                                                               bool value )
11.97.2.72 SetConfigurationBitBlu_LedToggleFast() void SetConfigurationBitBlu_LedToggleFast (
                                                                              bool value )
\textbf{11.97.2.73} \quad \textbf{SetConfigurationBitBlu\_LedToggleSlow()} \quad \texttt{void SetConfigurationBitBlu\_LedToggleSlow} \quad \textbf{(}
                                                                               bool value )
11.97.2.74 SetConfigurationBitCC_Gen() void SetConfigurationBitCC_Gen (
                                                                               bool value )
11.97.2.75 SetConfigurationBitCV_Gen() void SetConfigurationBitCV_Gen (
                                                                               bool value )
11.97.2.76 SetConfigurationBitRC_Gen() void SetConfigurationBitRC_Gen (
                                                                              bool value )
11.97.2.77 SetConfigurationBitRed_Led() void SetConfigurationBitRed_Led (
                                                                               bool value )
11.97.2.78 SetConfigurationBitRed LedSaturation() void SetConfigurationBitRed_LedSaturation (
                                                                               bool value )
11.97.2.79 SetConfigurationBitRed_LedToggleFast() void SetConfigurationBitRed_LedToggleFast (
                                                                              bool value )
\textbf{11.97.2.80} \quad \textbf{SetConfigurationBitRed\_LedToggleSlow()} \quad \texttt{void SetConfigurationBitRed\_LedToggleSlow} \quad \textbf{(} \quad \textbf{(} \quad \textbf{)} \quad \textbf{(} \quad \textbf{)} \quad \textbf{(} 
                                                                               bool value )
```

```
11.97.2.81 SetConfigurationBitRelais() void SetConfigurationBitRelais (
              bool value )
11.97.2.82 SetConfigurationBitRV_Gen() void SetConfigurationBitRV_Gen (
              bool value )
\textbf{11.97.2.83} \quad \textbf{SetConfigurationBitStream()} \quad \texttt{void SetConfigurationBitStream} \ \ \textbf{(}
              bool value )
11.97.2.84 SetConfigurationBitSupply() void SetConfigurationBitSupply (
              bool value )
11.97.2.85 SetCrossTalkOffset() void SetCrossTalkOffset (
              int CrossTalk )
11.97.2.86 SetCrossTalkOptimum() void SetCrossTalkOptimum (
              int cxOptimum )
11.97.2.87 SetDigout() void SetDigout (
              uint16_t index,
              bool enable )
11.97.2.88 SetDisplayText() void SetDisplayText (
              int index,
              String^{\wedge} displayText)
11.97.2.89 SetDownsampleFactor() void SetDownsampleFactor (
              int index,
              int downsample_factor )
```

```
11.97.2.90 SetFilter() void SetFilter (
             int filter )
11.97.2.91 SetFilterCoeffs() void SetFilterCoeffs (
             int index,
              array< int >^{\land} coeffs )
11.97.2.92 SetlClamp() void SetIClamp (
             int iClamp )
11.97.2.93 SetICOffset() void SetICOffset (
             int ICOffset )
11.97.2.94 SetlGain() void SetIGain (
             int iGain )
11.97.2.95 SetNoFilterCoeffs() void SetNoFilterCoeffs (
             int index )
11.97.2.96 SetPGain() void SetPGain (
             int pGain )
\textbf{11.97.2.97} \quad \textbf{SetRecordingNumber()} \quad \texttt{void SetRecordingNumber ()}
              unsigned int recordingNumber )
11.97.2.98 SetScreen() void SetScreen (
            int screen )
```


```
\textbf{11.97.2.99} \quad \textbf{SetSimulation()} \quad \texttt{void SetSimulation ()}
             unsigned int enable )
11.97.2.100 SetUClamp() void SetUClamp (
             int uClamp )
11.97.2.101 SetUCOffset() void SetUCOffset (
             int UCOffset )
11.97.2.102 SetUVOffset() void SetUVOffset (
             int UVOffset )
11.97.2.103 SetXGain() void SetXGain (
             int xGain )
11.97.2.104 StopTable() [1/2] void StopTable ( )
11.97.2.105 StopTable() [2/2] void StopTable (
             int timeout )
11.97.2.106 Table_Wait() void Table_Wait (
             unsigned int tableWait )
11.97.2.107 TableDefBegin() void TableDefBegin ( )
11.97.2.108 TableDefEnd() void TableDefEnd ( )
```

11.97.2.109 UpdateDisplay() void UpdateDisplay ()

11.98 CRoboDeviceNet Class Reference

CRoboDeviceNet is the base class for all Robo platform based devices

Inheritance diagram for CRoboDeviceNet:

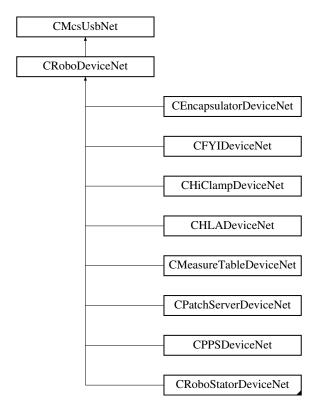

Classes

class RoboMainLowLevelCommands

Public Member Functions

- CRoboDeviceNet (void)
- ∼CRoboDeviceNet (void)
- void SetInMovement ()

Low level command, sets the internal state to "In Movement"

bool GetInMovement ()

Low level command, gets the internal state "In Movement"

uint32_t GetMovementError ()

Low level command, gets the error of the last movement end

- void FindReference (unsigned char busaddress, char axes)
- void FindReference (unsigned char busaddress, char axes, int timeout)

Searches the reference position of the motor

- void MoveAbs (unsigned char busaddress, char axes, int x, int y)
- void MoveAbs (unsigned char busaddress, char axes, int x, int y, int timeout)

Moves the motor to the new absolute position

- void StopMovement (unsigned char busaddress, char axes)
- void StopMovement (unsigned char busaddress, char axes, int timeout)

Stops the current movement

- · void SetCurrentAndAir (unsigned char busaddress, char axes, unsigned short onoff)
- void SetCurrentAndAir (unsigned char busaddress, char axes, unsigned short onoff, int timeout)
- void CancelPoolLoop ()
- void CancelPoolLoopAndStopMovement ()
- void GetCurrentPosition (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out]int%
 x, [System::Runtime::InteropServices::Out]int% y)

Gets the current position of motors

- · void SetAirValve (unsigned int onoff)
- unsigned int GetAirValve ()
- unsigned int GetVoltageValves ()
- unsigned int GetVoltageRs485A ()
- unsigned int GetVoltageRs485B ()
- · unsigned int GetVoltageAirvalve ()
- unsigned int GetCurrentAirvalve ()
- unsigned int GetVoltage12V ()
- unsigned int GetAirpressure ()
- unsigned int GetVoltage5V ()
- unsigned int GetErrorVoltageValves ()
- unsigned int GetErrorVoltageRs485A ()
- unsigned int GetErrorVoltageRs485B ()
- · unsigned int GetErrorVoltageAirvalve ()
- unsigned int GetErrorCurrentAirvalve ()
- unsigned int GetErrorVoltage12V ()
- unsigned int GetErrorAirpressure ()
- unsigned int GetErrorVoltage5V ()
- void SetVoltageValvesLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageRs485ALimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageRs485BLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageAirvalveLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetCurrentAirvalveLimit (unsigned int lowercurrent, unsigned int uppercurrent)
- void SetVoltage12VLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetAirpressureLimit (unsigned int lowerpressure, unsigned int upperpressure)
- void SetVoltage5VLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void GetVoltageValvesLimit ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System←::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltageRs485ALimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltageRs485BLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltageAirvalveLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetCurrentAirvalveLimit ([System::Runtime::InteropServices::Out]unsigned int% lowercurrent, [System::Runtime::InteropServices::Out]unsigned int% uppercurrent)
- void GetVoltage12VLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::
 Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetAirpressureLimit ([System::Runtime::InteropServices::Out]unsigned int% lowerpressure, [System ← ::Runtime::InteropServices::Out]unsigned int% upperpressure)
- void GetVoltage5VLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::
 Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void SetMinPressure (int pressure)
- int GetMinPressure ()

Static Public Attributes

```
    static const uint32 t RoboError Base = (0xA0110000L)

    static const uint32 t RoboError UnknownCommand = ( (0xA0110000L) )

static const uint32_t RoboError_Timeout = ( (0xA0110000L) | 0x0001 )

    static const uint32 t RoboError Pressure = ( (0xA0110000L) | 0x0002 )

    static const uint32 t RoboError RangeExceeded = ( (0xA0110000L) | 0x0003 )

    static const uint32_t RoboError_CommunicationTimeout = ( (0xA0110000L) | 0x0004 )

    static const uint32_t RoboError_AnotherMaster = ( (0xA0110000L) | 0x0005 )

    static const uint32 t RoboError FindReferenceMethod = ( (0xA0110000L) | 0x0006 )

    static const uint32 t RoboError NoSpeedOrAcceleration = ( (0xA0110000L) | 0x0007 )

    static const uint32 t RoboError NoEndSwitch = ( (0xA0110000L) | 0x0008 )

    static const uint32_t RoboError_CannotEscapeEndSwitch = ( (0xA0110000L) | 0x0009 )

    static const uint32 t RoboError CommandAlreadyInProgress = ( (0xA0110000L) | 0x000A )

    static const uint32_t RoboError_NoReference = ( (0xA0110000L) | 0x000B )

    static const uint32 t RoboError OverPressure = ( (0xA0110000L) | 0x000C )

    static const uint32 t RoboError Phase0OutOfRange = ( (0xA0110000L) | 0x000D )

    static const uint32_t RoboError_PeristalticTimeout = ( (0xA0110000L) | 0x000E )

    static const uint32 t RoboError GilsonTimeout = ( (0xA0110000L) | 0x000F )

    static const uint32_t RoboError_GilsonWrondID = ( (0xA0110000L) | 0x0010 )

    static const uint32 t RoboError GilsonCommandPending = ((0xA0110000L) | 0x0011)

    static const uint32 t RoboError ParameterNotAllowed = ( (0xA0110000L) | 0x0012 )

    static const uint32_t RoboError_StateChangeNotPossible = ( (0xA0110000L) | 0x0013 )

    static const uint32 t RoboError CommandNotPossible = ( (0xA0110000L) | 0x0014 )

static const uint32_t RoboError_DacqNotReady = ( (0xA0110000L) | 0x0015 )
static const uint32_t RoboError_NoMoreData = ( (0xA0110000L) | 0x0016 )

    static const uint32 t RoboError McsBus UnknownCommand = ( (0xA0110000L) | 0x003F)

    static const uint32 t RoboError DLLMovementTimeout = ( (0xA0110000L) | 0x1001)

    static const uint32 t RoboError PollLoopCanceled = ( (0xA0110000L) | 0x1002)

    static const uint32_t RoboError_PollLoopCanceledAndStopMovement = ( (0xA0110000L) | 0x1003)

• static const byte McsBus_XY = 1
     McsBus address for the xy-plane

    static const byte McsBus ZI = 2

     McsBus address for the z and i axes
• static const byte Axis X = 0
     Axis number of x for axis argument

    static const byte Axis_Y = 1

     Axis number of y for axis argument
• static const byte Axis Z = 0
     Axis number of z for axis argument

    static const byte Axis I = 1

     Axis number of i for axis argument

    static const char Axes X = 1

     Bit pattern for x axis for axes argument

 static const char Axes Y = 2

     Bit pattern for y axis for axes argument

    static const char Axes Z = 1

     Bit pattern for z axis for axes argument
static const char Axes_I = 2
     Bit pattern for i axis for axes argument
```

Properties

- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet^{\(\Lambda\)} McsBus_MotorControl [get]
- RoboMainLowLevelCommands RoboMainLowLevelCommand [get]

Events

RoboStatusEventDelegate[^] RoboStatusEvent

Additional Inherited Members

11.98.1 Detailed Description

CRoboDeviceNet is the base class for all Robo platform based devices

11.98.2 Constructor & Destructor Documentation

```
11.98.2.1 CRoboDeviceNet() CRoboDeviceNet (
void )

11.98.2.2 ~CRoboDeviceNet() ~CRoboDeviceNet (
void )
```

11.98.3 Member Function Documentation

```
11.98.3.1 CancelPoolLoop() void CancelPoolLoop ( )
```

 $\textbf{11.98.3.2} \quad \textbf{CancelPoolLoopAndStopMovement()} \quad \texttt{void CancelPoolLoopAndStopMovement ()} \\$

```
11.98.3.3 FindReference() [1/2] void FindReference (
unsigned char busaddress,
char axes)
```

Searches the reference position of the motor

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)

```
11.98.3.5 GetAirpressure() unsigned int GetAirpressure ( )
```

```
11.98.3.6 GetAirpressureLimit() void GetAirpressureLimit (

[System::Runtime::InteropServices::Out] unsigned int% lowerpressure,

[System::Runtime::InteropServices::Out] unsigned int% upperpressure)
```

```
11.98.3.7 GetAirValve() unsigned int GetAirValve ( )
```

```
11.98.3.8 GetCurrentAirvalve() unsigned int GetCurrentAirvalve ( )
```

Gets the current position of motors

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
X	Current position of first axis if pattern in axes is set
У	Current position of second axis if pattern in axes is set

11.98.3.11	GetErrorAirpressure() unsigned int GetErrorAirpressure ()
11.98.3.12	GetErrorCurrentAirvalve() unsigned int GetErrorCurrentAirvalve ()
11.98.3.13	GetErrorVoltage12V() unsigned int GetErrorVoltage12V ()
11.98.3.14	<pre>GetErrorVoltage5V() unsigned int GetErrorVoltage5V ()</pre>
11.98.3.15	<pre>GetErrorVoltageAirvalve() unsigned int GetErrorVoltageAirvalve ()</pre>
11.98.3.16	GetErrorVoltageRs485A() unsigned int GetErrorVoltageRs485A ()
11.98.3.17	<pre>GetErrorVoltageRs485B() unsigned int GetErrorVoltageRs485B ()</pre>
11.98.3.18	GetErrorVoltageValves() unsigned int GetErrorVoltageValves ()
	GetInMovement() bool GetInMovement () ommand, gets the internal state "In Movement"
11.98.3.20	<pre>GetMinPressure() int GetMinPressure ()</pre>

```
11.98.3.21 GetMovementError() uint32_t GetMovementError ( )
Low level command, gets the error of the last movement end
11.98.3.22 GetVoltage12V() unsigned int GetVoltage12V ( )
11.98.3.23 GetVoltage12VLimit() void GetVoltage12VLimit (
             [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
             [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
11.98.3.24 GetVoltage5V() unsigned int GetVoltage5V ()
11.98.3.25 GetVoltage5VLimit() void GetVoltage5VLimit (
             [{\tt System::Runtime::InteropServices::Out}] \ {\tt unsigned int \$} \ {\tt lowervoltage,}
             [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
11.98.3.26 GetVoltageAirvalve() unsigned int GetVoltageAirvalve ()
11.98.3.27 GetVoltageAirvalveLimit() void GetVoltageAirvalveLimit (
             [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
             [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
11.98.3.28 GetVoltageRs485A() unsigned int GetVoltageRs485A ()
11.98.3.29 GetVoltageRs485ALimit() void GetVoltageRs485ALimit (
             [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
             [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
```

```
\textbf{11.98.3.30} \quad \textbf{GetVoltageRs485B()} \quad \texttt{unsigned int GetVoltageRs485B ()}
11.98.3.31 GetVoltageRs485BLimit() void GetVoltageRs485BLimit (
              [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
              [{\tt System::Runtime::InteropServices::Out}] \  \, unsigned \  \, int {\tt \textit{uppervoltage}} \  \, )
11.98.3.32 GetVoltageValves() unsigned int GetVoltageValves ( )
11.98.3.33 GetVoltageValvesLimit() void GetVoltageValvesLimit (
              [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
              [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
11.98.3.34 MoveAbs() [1/2] void MoveAbs (
              unsigned char busaddress,
              char axes,
              int x,
              int y)
11.98.3.35 MoveAbs() [2/2] void MoveAbs (
              unsigned char busaddress,
              char axes,
              int x,
              int y,
              int timeout )
```

Moves the motor to the new absolute position

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
X	Position of first axis, if pattern in axes is set
У	Position of second axis if pattern in axes is set
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)

11.98.3.36 SetAirpressureLimit() void SetAirpressureLimit (

```
unsigned int lowerpressure,
             unsigned int upperpressure )
11.98.3.37 SetAirValve() void SetAirValve (
             unsigned int onoff )
11.98.3.38 SetCurrentAirvalveLimit() void SetCurrentAirvalveLimit (
             unsigned int lowercurrent,
             unsigned int uppercurrent )
11.98.3.39 SetCurrentAndAir() [1/2] void SetCurrentAndAir (
             unsigned char busaddress,
             char axes,
             unsigned short onoff )
11.98.3.40 SetCurrentAndAir() [2/2] void SetCurrentAndAir (
             unsigned char busaddress,
             char axes,
             unsigned short onoff,
             int timeout )
11.98.3.41 SetInMovement() void SetInMovement ()
Low level command, sets the internal state to "In Movement"
11.98.3.42 SetMinPressure() void SetMinPressure (
             int pressure )
11.98.3.43 SetVoltage12VLimit() void SetVoltage12VLimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
```

```
11.98.3.44 SetVoltage5VLimit() void SetVoltage5VLimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
11.98.3.45 SetVoltageAirvalveLimit() void SetVoltageAirvalveLimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
11.98.3.46 SetVoltageRs485ALimit() void SetVoltageRs485ALimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
11.98.3.47 SetVoltageRs485BLimit() void SetVoltageRs485BLimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
11.98.3.48 SetVoltageValvesLimit() void SetVoltageValvesLimit (
             unsigned int lowervoltage,
             unsigned int uppervoltage )
11.98.3.49 StopMovement() [1/2] void StopMovement (
             unsigned char busaddress,
             char axes )
11.98.3.50 StopMovement() [2/2] void StopMovement (
             unsigned char busaddress,
             char axes,
             int timeout )
```

Stops the current movement

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)

11.98.4 Member Data Documentation

11.98.4.1 Axes_I const char Axes_I = 2 [static]

Bit pattern for i axis for axes argument

11.98.4.2 Axes_X const char Axes_X = 1 [static]

Bit pattern for x axis for axes argument

11.98.4.3 Axes_Y const char Axes_Y = 2 [static]

Bit pattern for y axis for axes argument

11.98.4.4 Axes_Z const char Axes_Z = 1 [static]

Bit pattern for z axis for axes argument

11.98.4.5 Axis_I const byte Axis_I = 1 [static]

Axis number of i for axis argument

11.98.4.6 Axis_X const byte $Axis_X = 0$ [static]

Axis number of x for axis argument

11.98.4.7 Axis_Y const byte Axis_Y = 1 [static]

Axis number of y for axis argument

```
11.98.4.8 Axis_Z const byte Axis_Z = 0 [static]
```

Axis number of z for axis argument

```
11.98.4.9 McsBus_XY const byte McsBus_XY = 1 [static]
```

McsBus address for the xy-plane

```
11.98.4.10 McsBus_ZI const byte McsBus_ZI = 2 [static]
```

McsBus address for the z and i axes

```
11.98.4.11 RoboError_AnotherMaster const uint32_t RoboError_AnotherMaster = ( (0xA0110000L) | 0x0005 ) [static]
```

```
11.98.4.12 RoboError_Base const uint32_t RoboError_Base = (0xA0110000L) [static]
```

11.98.4.13 RoboError_CannotEscapeEndSwitch const uint32_t RoboError_CannotEscapeEndSwitch = ((0xA0110000L) | 0x0009) [static]

11.98.4.14 RoboError_CommandAlreadyInProgress const uint32_t RoboError_CommandAlreadyIn↔ Progress = ((0xA0110000L) | 0x000A) [static]

11.98.4.15 RoboError_CommandNotPossible const uint32_t RoboError_CommandNotPossible = ((0xA0110000L) | 0x0014) [static]

11.98.4.16 RoboError_CommunicationTimeout const uint32_t RoboError_CommunicationTimeout = ((0xA0110000L) | 0x0004) [static]

```
11.98.4.17 RoboError_DacqNotReady const uint32_t RoboError_DacqNotReady = ( (0xA0110000L) |
0x0015 ) [static]
11.98.4.18 RoboError_DLLMovementTimeout const uint32_t RoboError_DLLMovementTimeout = (
(0xA0110000L) | 0x1001) [static]
11.98.4.19 RoboError_FindReferenceMethod const uint32_t RoboError_FindReferenceMethod = (
(0xA0110000L) | 0x0006 ) [static]
11.98.4.20 RoboError_GilsonCommandPending const uint32_t RoboError_GilsonCommandPending = (
(0xA0110000L) | 0x0011 ) [static]
11.98.4.21 RoboError_GilsonTimeout const uint32_t RoboError_GilsonTimeout = ( (0xA0110000L) |
0x000F ) [static]
11.98.4.22 RoboError_GilsonWrondID const uint32_t RoboError_GilsonWrondID = ( (0xA0110000L) |
0x0010 ) [static]
11.98.4.23 RoboError_McsBus_UnknownCommand const uint32_t RoboError_McsBus_UnknownCommand
= ((0xA0110000L) | 0x003F) [static]
11.98.4.24 RoboError NoEndSwitch const uint32_t RoboError_NoEndSwitch = ( (0xA0110000L) |
0x0008 ) [static]
11.98.4.25 RoboError_NoMoreData const uint32_t RoboError_NoMoreData = ( (0xA0110000L) | 0x0016
) [static]
11.98.4.26 RoboError_NoReference const uint32_t RoboError_NoReference = ( (0xA0110000L) |
0x000B ) [static]
```

```
11.98.4.27 RoboError_NoSpeedOrAcceleration const uint32_t RoboError_NoSpeedOrAcceleration = (
  (0xA0110000L) | 0x0007 ) [static]
11.98.4.28 RoboError_OverPressure const uint32_t RoboError_OverPressure = ( (0xA0110000L) |
0x000C ) [static]
11.98.4.29 RoboError_ParameterNotAllowed const uint32_t RoboError_ParameterNotAllowed = (
  (0xA0110000L) | 0x0012 ) [static]
11.98.4.30 RoboError_PeristalticTimeout const uint32_t RoboError_PeristalticTimeout = ( (0x←
A0110000L) | 0x000E ) [static]
\textbf{11.98.4.31} \quad \textbf{RoboError\_Phase0OutOfRange} \quad \texttt{const uint} \\ \textbf{32\_t RoboError\_Phase0OutOfRange} = \textbf{((0x} \leftarrow \texttt{0x}) \\ \textbf{(0x} \leftarrow \texttt{0x}) 
A0110000L) | 0x000D ) [static]
11.98.4.32 RoboError_PollLoopCanceled const uint32_t RoboError_PollLoopCanceled = ( (0x←
A0110000L) | 0x1002) [static]
11.98.4.33 RoboError_PollLoopCanceledAndStopMovement const uint32_t RoboError_PollLoop←
{\tt CanceledAndStopMovement = ( (0xA0110000L) | 0x1003) [static]}
11.98.4.34 RoboError Pressure const uint32_t RoboError_Pressure = ( (0xA0110000L) | 0x0002 )
 [static]
11.98.4.35 RoboError_RangeExceeded const uint32_t RoboError_RangeExceeded = ( (0xA0110000L) |
0x0003 ) [static]
11.98.4.36 RoboError_StateChangeNotPossible const uint32_t RoboError_StateChangeNotPossible =
  ( (0xA0110000L) | 0x0013 ) [static]
```

```
11.98.4.37 RoboError_Timeout const uint32_t RoboError_Timeout = ( (0xA0110000L) | 0x0001) [static]
```

```
11.98.4.38 RoboError_UnknownCommand const uint32_t RoboError_UnknownCommand = ( (0x \leftarrow A0110000L) ) [static]
```

11.98.5 Property Documentation

```
11.98.5.1 McsBus CMcsBusNet^ McsBus [get]
```

11.98.5.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]


11.98.5.3 RoboMainLowLevelCommand RoboMainLowLevelCommands^ RoboMainLowLevelCommand [get]

11.98.6 Event Documentation

11.98.6.1 RoboStatusEvent RoboStatusEventDelegate^ RoboStatusEvent

11.99 CRoboFluidDeviceNet Class Reference

Inheritance diagram for CRoboFluidDeviceNet:

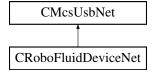

Public Member Functions

- CRoboFluidDeviceNet (void)
- ∼CRoboFluidDeviceNet (void)
- void SetValve (int value)

Open or Close valves.

void SetSingleValve (int valve, bool onoff)

Opens or Closes a valve.

• int GetValve ()

Query the state of the values.

bool GetSingleValve (int valve)

Query the state of a valve.

- void CloseAllValves ()
- void PumpOn (int index, short speed)
- void SetPumpSpeed (int index, short speed)
- void PumpOff (int index)
- short GetPumpSpeed (int index)
- bool IsPumpMotorOn (int index)

Protected Attributes

- CRoboFluidDevice * m_pRoboFluidDevice
- CMcsBus_MotorControlNet ^ m_pMcsBus_MotorControlNet

Properties

CMcsBus MotorControlNet[^] McsBus MotorControl [get]

Additional Inherited Members

11.99.1 Constructor & Destructor Documentation

```
11.99.1.1 CRoboFluidDeviceNet() CRoboFluidDeviceNet (
void )
```

```
11.99.1.2 \sim CRoboFluidDeviceNet() \sim CRoboFluidDeviceNet ( void )
```

11.99.2 Member Function Documentation

```
11.99.2.1 CloseAllValves() void CloseAllValves ( )
11.99.2.2 GetPumpSpeed() short GetPumpSpeed (
              int index )
11.99.2.3 GetSingleValve() bool GetSingleValve (
              int valve )
Query the state of a valve.
Parameters
 valve
        number of valve /*!
Returns
     state of the valve
11.99.2.4 GetValve() int GetValve ()
Query the state of the values.
Returns
     the current state of the valves as a bit pattern.
11.99.2.5 IsPumpMotorOn() bool IsPumpMotorOn (
             int index )
11.99.2.6 PumpOff() void PumpOff (
             int index )
```

11.99.2.7 PumpOn() void PumpOn (
 int index,
 short speed)

Opens or Closes a valve.

Parameters

valve	number of valve to be changed /*!
onoff	open or close the valve

```
11.99.2.10 SetValve() void SetValve ( int value )
```

Open or Close valves.

Parameters

Ī	value	bit pattern of valves which should be open.
---	-------	---

11.99.3 Member Data Documentation

11.99.3.1 m_pMcsBus_MotorControlNet CMcsBus_MotorControlNet ^ m_pMcsBus_MotorControlNet [protected]

11.99.3.2 m_pRoboFluidDevice CRoboFluidDevice* m_pRoboFluidDevice [protected]


11.99.4 Property Documentation

11.99.4.1 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.100 CRobolnjectDeviceNet Class Reference

CRobolnjectDeviceNet is the to control the MCS Robolnject device

Inheritance diagram for CRobolnjectDeviceNet:

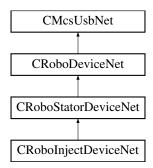

Public Member Functions

· CRobolnjectDeviceNet (void)

Additional Inherited Members

11.100.1 Detailed Description

CRobolnjectDeviceNet is the to control the MCS Robolnject device


11.100.2 Constructor & Destructor Documentation

```
11.100.2.1 CRobolnjectDeviceNet() CRobolnjectDeviceNet (
void )
```

11.101 CRoboocyte2DeviceNet Class Reference

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

Inheritance diagram for CRoboocyte2DeviceNet:

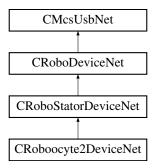

Public Member Functions

- CRoboocyte2DeviceNet (void)
- void SetAxisLED (bool onoff)
- bool GetAxisLED ()
- CRoboDacqNet ^ GetRoboDacq ()
- CRoboFluidDeviceNet ^ GetRoboFluidDevice ()
- CGilsonDeviceNet ^ GetGilsonDevice ()
- CMcsBus_ExtensionNet ^ GetMcsBus_Extension ()

Additional Inherited Members

11.101.1 Detailed Description

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

11.101.2 Constructor & Destructor Documentation

```
11.101.2.1 CRoboocyte2DeviceNet() CRoboocyte2DeviceNet (
void )
```

11.101.3 Member Function Documentation

```
11.101.3.1 GetAxisLED() bool GetAxisLED ( )
```

```
11.101.3.2 GetGilsonDevice() CGilsonDeviceNet ^ GetGilsonDevice ( )
```

```
11.101.3.3 GetMcsBus_Extension() CMcsBus_ExtensionNet ^ GetMcsBus_Extension ( )
```


11.101.3.4 GetRoboDacq() CRoboDacqNet ^ GetRoboDacq ()

11.101.3.5 GetRoboFluidDevice() CRoboFluidDeviceNet ^ GetRoboFluidDevice ()

```
11.101.3.6 SetAxisLED() void SetAxisLED (
bool onoff)
```

11.102 CRoboStatorDeviceNet Class Reference

Inheritance diagram for CRoboStatorDeviceNet:

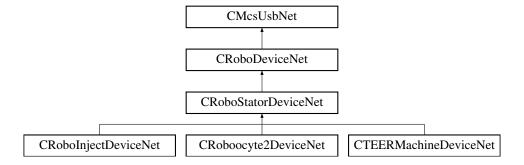

Classes

class RoboMainStatorLowLevelCommands

Public Member Functions

- CRoboStatorDeviceNet (void)
- void FindReferenceXY ()
- void FindReferenceXY (int timeout)
- void FindReferenceZ ()
- void FindReferenceZ (int timeout)
- void FindReferencel ()
- void FindReferencel (int timeout)
- unsigned char HasRefXY ()
- unsigned char HasRefZ ()
- unsigned char HasRefl ()
- void MoveAbsXY (int x, int y)
- void MoveAbsXY (int x, int y, int timeout)
- void MoveAbsZ (int z)
- void MoveAbsZ (int z, int timeout)
- void MoveAbsI (int i)
- void MoveAbsI (int i, int timeout)
- void StopMovementXY ()
- void StopMovementXY (int timeout)
- void StopMovementZ ()
- void StopMovementZ (int timeout)
- void StopMovementI ()
- void StopMovementI (int timeout)
- void SetCurrentAndAirXY (unsigned short onoff)

- void SetCurrentAndAirXY (unsigned short onoff, int timeout)
- void GetCurrentPositionXY ([System::Runtime::InteropServices::Out]int% x, [System::Runtime::Interop←
 Services::Out]int% y)
- void GetCurrentPositionZ ([System::Runtime::InteropServices::Out]int% z)
- void GetCurrentPositionI ([System::Runtime::InteropServices::Out]int% i)
- void SetVelocityXY (int v)
- void SetVelocityZ (int v)
- void SetVelocityI (int v)
- void SetSpeedXY (int v)
- void SetSpeedZ (int v)
- void SetSpeedI (int v)
- void SetSpeedNativeXY (int v)
- void SetSpeedNativeZ (int v)
- void SetSpeedNativel (int v)
- void SetAccelerationXY (int a)
- void SetAccelerationZ (int a)
- void SetAccelerationI (int a)
- void SetAccelerationNativeXY (int a)
- void SetAccelerationNativeZ (int a)
- void SetAccelerationNativel (int a)

Properties

• RoboMainStatorLowLevelCommands RoboMainStatorLowLevelCommand [get]

Additional Inherited Members

11.102.1 Constructor & Destructor Documentation

```
11.102.1.1 CRoboStatorDeviceNet() CRoboStatorDeviceNet (
```

11.102.2 Member Function Documentation

```
11.102.2.1 FindReferencel() [1/2] void FindReferenceI ( )
```

```
11.102.2.2 FindReferencel() [2/2] void FindReferenceI ( int timeout )
```

```
11.102.2.3 FindReferenceXY() [1/2] void FindReferenceXY ( )
11.102.2.4 FindReferenceXY() [2/2] void FindReferenceXY (
             int timeout )
11.102.2.5 FindReferenceZ() [1/2] void FindReferenceZ ( )
11.102.2.6 FindReferenceZ() [2/2] void FindReferenceZ (
             int timeout )
11.102.2.7 GetCurrentPositionI() void GetCurrentPositionI (
             [System::Runtime::InteropServices::Out] int% i )
11.102.2.8 GetCurrentPositionXY() void GetCurrentPositionXY (
             [System::Runtime::InteropServices::Out] int% x,
             [System::Runtime::InteropServices::Out] int% y )
11.102.2.9 GetCurrentPositionZ() void GetCurrentPositionZ (
             [System::Runtime::InteropServices::Out] int% z )
11.102.2.10 HasRefl() unsigned char HasRefI ()
11.102.2.11 HasRefXY() unsigned char HasRefXY ( )
11.102.2.12 HasRefZ() unsigned char HasRefZ ( )
```

```
11.102.2.13 MoveAbsl() [1/2] void MoveAbsI (
             int i)
11.102.2.14 MoveAbsI() [2/2] void MoveAbsI (
             int i,
             int timeout )
11.102.2.15 MoveAbsXY() [1/2] void MoveAbsXY (
             int x,
             int y )
11.102.2.16 MoveAbsXY() [2/2] void MoveAbsXY (
             int x,
              int y,
             int timeout )
11.102.2.17 MoveAbsZ() [1/2] void MoveAbsZ (
             int z )
11.102.2.18 MoveAbsZ() [2/2] void MoveAbsZ (
             int z,
             int timeout )
11.102.2.19 SetAccelerationI() void SetAccelerationI (
             int a )
\textbf{11.102.2.20} \quad \textbf{SetAccelerationNativeI()} \quad \texttt{void SetAccelerationNativeI} \quad (
             int a)
11.102.2.21 SetAccelerationNativeXY() void SetAccelerationNativeXY (
             int a)
```

```
\textbf{11.102.2.22} \quad \textbf{SetAccelerationNativeZ()} \quad \texttt{void SetAccelerationNativeZ} \quad (
               int a)
11.102.2.23 SetAccelerationXY() void SetAccelerationXY (
              int a )
11.102.2.24 SetAccelerationZ() void SetAccelerationZ (
              int a)
11.102.2.25 SetCurrentAndAirXY() [1/2] void SetCurrentAndAirXY (
               unsigned short onoff )
11.102.2.26 SetCurrentAndAirXY() [2/2] void SetCurrentAndAirXY (
               unsigned short onoff,
               int timeout )
11.102.2.27 SetSpeedI() void SetSpeedI (
              int v)
11.102.2.28 SetSpeedNativel() void SetSpeedNativeI (
               int v)
11.102.2.29 SetSpeedNativeXY() void SetSpeedNativeXY (
              int v )
\textbf{11.102.2.30} \quad \textbf{SetSpeedNativeZ()} \quad \texttt{void SetSpeedNativeZ} \quad \textbf{(}
               int v)
```

```
\textbf{11.102.2.31} \quad \textbf{SetSpeedXY()} \quad \texttt{void SetSpeedXY} \ \ (
             int v )
11.102.2.32 SetSpeedZ() void SetSpeedZ (
             int v)
11.102.2.33 SetVelocityI() void SetVelocityI (
             int v)
11.102.2.34 SetVelocityXY() void SetVelocityXY (
             int v)
11.102.2.35 SetVelocityZ() void SetVelocityZ (
             int v)
11.102.2.36 StopMovementI() [1/2] void StopMovementI ( )
11.102.2.37 StopMovementI() [2/2] void StopMovementI (
             int timeout )
11.102.2.38 StopMovementXY() [1/2] void StopMovementXY ( )
11.102.2.39 StopMovementXY() [2/2] void StopMovementXY (
             int timeout )
11.102.2.40 StopMovementZ() [1/2] void StopMovementZ ( )
```


11.102.2.41 StopMovementZ() [2/2] void StopMovementZ (int timeout)

11.102.3 Property Documentation

11.102.3.1 RoboMainStatorLowLevelCommand RoboMainStatorLowLevelCommands ↑ RoboMainStator← LowLevelCommand [get]

11.103 CSafeISDeviceNet Class Reference

Inheritance diagram for CSafeISDeviceNet:

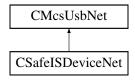

Public Member Functions

CSafeISDeviceNet (void)

Initializes a new instance of the CSafeISDeviceNet class.

∼CSafeISDeviceNet (void)

Releases unmanaged resources and performs other cleanup operations before the CSafelSDeviceNet is reclaimed by garbage collection.

void SetSwitches (unsigned short switches)

Sets the switches for all electrodes on the device. Do not use during measurement

void SetAdcChannels (unsigned char channels)

Sets the ADC channels you want to be sampled

void SetAdcSamplePos (array< unsigned short >^ positions)

Sets the sample position of the ADC.

• void SetDacMode (unsigned char mode)

Sets the DAC mode.

void SetDacPulseform (array< short >^ pulseform)

Sets the DAC pulseform.

• void SetDacPeriode (unsigned int periode)

Sets the DAC periode.

Properties

CRoboDeviceNet[^] RoboDevice [get]

Gets the CRoboDeviceNet. Use this to control the syringe.

• CFluidControlDeviceNet^ FluidControlDevice [get]

Gets the CFluidControlDeviceNet. Use this to control the valves. Only SetSingleValve is implemented for CSafeISDeviceNet.

CMcsUsbDacqNet[^] DacqDevice [get]

Gets the CMcsUsbDacqNet. Use this to control the data aquisition.

Additional Inherited Members

11.103.1 Detailed Description

11.103.2 Constructor & Destructor Documentation

```
11.103.2.1 CSafeISDeviceNet() CSafeISDeviceNet (
```

Initializes a new instance of the CSafeISDeviceNet class.

```
11.103.2.2 \simCSafeISDeviceNet() \simCSafeISDeviceNet ( void )
```

Releases unmanaged resources and performs other cleanup operations before the CSafelSDeviceNet is reclaimed by garbage collection.

11.103.3 Member Function Documentation

```
11.103.3.1 SetAdcChannels() void SetAdcChannels (
unsigned char channels )
```

Sets the ADC channels you want to be sampled

Parameters

channels The bitmap of the 8 channels. Set bit to 1 for the channels you want measure

Sets the sample position of the ADC.

Parameters

positions The positions in units of 0.1µs.

```
11.103.3.3 SetDacMode() void SetDacMode ( unsigned char mode )
```

Sets the DAC mode.

Parameters

```
mode | The mode: 0 = Impedance; 1 = Amperometry
```

Sets the DAC periode.

Parameters

periode The periode in units of 10μs.

Sets the DAC pulseform.

Parameters

```
pulseform The pulseform.
```

```
11.103.3.6 SetSwitches() void SetSwitches (
          unsigned short switches )
```

Sets the switches for all electrodes on the device. Do not use during measurement

Parameters

switches The switches: See Schematics for the meaning

11.103.4 Property Documentation

```
11.103.4.1 DacqDevice CMcsUsbDacqNet^ DacqDevice [get]
```

Gets the CMcsUsbDacqNet. Use this to control the data aquisition.

```
11.103.4.2 FluidControlDevice CFluidControlDeviceNet^ FluidControlDevice [get]
```


Gets the CFluidControlDeviceNet. Use this to control the valves. Only SetSingleValve is implemented for CSafeISDeviceNet.

```
11.103.4.3 RoboDevice CRoboDeviceNet^ RoboDevice [get]
```

Gets the CRoboDeviceNet. Use this to control the syringe.

11.104 CSCUDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CSCUDacqGroupChannelSelectionNet:

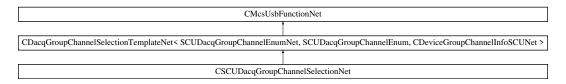

Public Member Functions

CSCUDacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members


11.104.1 Constructor & Destructor Documentation

11.104.1.1 CSCUDacqGroupChannelSelectionNet() CSCUDacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)

11.105 CSCUFunctionNet Class Reference

CSCUFunctionNet is the class to control the SCU device

Inheritance diagram for CSCUFunctionNet:

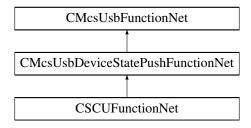

Public Member Functions

- delegate void OnGetAvailableHeadstages (uint32_t AvailableHeadstages)
- delegate void OnlsHeadstageAvailable (uint32_t Headstage, bool available)
- CSCUFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pSCUFunctionPointer
 — Container)

Initializes a new instance of the CSCUFunctionNet class.

- CSCUFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CSCUFunctionNet ()
- !CSCUFunctionNet ()
- uint32_t GetAvailableHeadstages ()

Gets a bitmap of available headstages.

bool IsInDacqLegacyMode ()

Is the SCU in legacy mode

void SetDacqLegacyMode (bool enable)

Enable the SCU legacy mode

• uint32_t GetMaxStimulusChannelsPerHeadstage ()

Gets the maximal number of stimulation channels a headstage can have.

uint32_t GetMaxNumberOfHeadstages ()

Gets the maximal number of headstages.

SCU_HeadstageIdEnumNet GetHeadstageID (uint32_t Headstage)

Gets the headstage fpga ID.

bool IsHeadstageAvailable (uint32_t Headstage)

Checks whether the given headstage is available.

void PowerHS (uint32 t Headstage, bool power)

Power the HS

• bool IsHSPowered (uint32_t Headstage)

Is the HS powered

bool HasHSPowerSwitch ()

Has SCU HS power switch

String \(^\) GetHeadstageSerialNumber (uint32 t Headstage)

Gets the serial number of a given headstage.

uint32_t GetHeadstageNumberOfAnalogChannels (uint32_t Headstage)

Gets the number of analog channels for a given headstage.

• uint32_t GetHeadstageNumberOfStimulationChannels (uint32_t Headstage)

Gets the number of stimulation channels for a given headstage.

uint32_t GetHeadstageGainInPermille (uint32_t Headstage)

Gets the gain factor in permille for a given headstage.

uint32 t GetHeadstageAdcRangeInMicroVolt (uint32 t Headstage)

Gets the ADC Range in uV for a given headstage.

uint32 t GetHeadstageAdcBits (uint32 t Headstage)

Gets the Number of ADC bits for a given headstage.

uint32_t GetHeadstageDacVoltageRangeInMilliVolt (uint32_t Headstage)

Gets the DAC Voltage Range in mV for a given headstage.

uint32 t GetHeadstageDacVoltageResolutionInMicroVolt (uint32 t Headstage)

Gets the DAC Voltage Resolution in uV for a given headstage.

uint32_t GetHeadstageDacCurrentRangeInMicroAmpere (uint32_t Headstage)

Gets the DAC Current Range in uA for a given headstage.

uint32_t GetHeadstageDacCurrentResolutionInNanoAmpere (uint32_t Headstage)

Gets the DAC Current Resolution in nA for a given headstage.

• uint32_t GetHeadstageDacBits (uint32_t Headstage)

Gets the Number of DAC bits for a given headstage.

uint32 t GetHeadstageSamplerate (uint32 t Headstage)

Gets the Samplerate of a given headstage.

bool HasGalvanicIsolation ()

Has galvanic isolated hardware

bool HasAnalogOut ()

Has AnalogOut hardware

void EnableAnalogOut (bool enable)

Enables AnalogOut globally

bool IsAnalogOutEnabled ()

Is AnalogOut enabled

void SetAnalogOutDACRange (AnalogOut_DAC_Range_EnumNet range)

Sets the analog out DAC range

AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ()

Gets the analog out DAC range

void SetAnalogOutADCRange (uint32_t range)

Sets the analog out ADC range

uint32_t GetAnalogOutADCRange ()

Gets the analog out ADC range

void AutomaticAnalogOut (bool automatic)

Sets automatic source channel selection

bool IsAutomaticAnalogOut ()

Is Automatic source channel selection selected

void SetAnalogOutChannels (uint32 t out channel, uint32 t source channel)

Set the source channel number for a certain output channel

uint32_t GetAnalogOutChannels (uint32_t out_channel)

Get the connected source channel number for a certain output channel

 void SetReferenceElectrodeSwitchState (uint32_t Headstage, ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)

Sets the position of the switch for the reference electrode

ReferenceElectrodeSwitchPositionEnumNet GetReferenceElectrodeSwitchState (uint32_t Headstage)

Gets the position of the switch for the reference electrode

void SetReferenceElectrodeMode (uint32_t Headstage, ReferenceElectrodeModeEnumNet NewValue)

Sets the mode for the reference electrode

• ReferenceElectrodeModeEnumNet GetReferenceElectrodeMode (uint32_t Headstage)

Gets the mode for the reference electrode

• CFilterPropertyNet ^ GetFilterProperty (SCUDacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)

Gets the filter property

array< CFilterPropertyNet[^]> [^] GetFilterProperties (SCUDacqGroupChannelEnumNet GroupID, int filter
 — Configurations_Length)

Gets multiple filter properties

Events

- OnGetAvailableHeadstages^ GetAvailableHeadstagesEvent [add, remove, raise]

 Event fires when the bitmap of available headstages has changed
- OnlsHeadstageAvailable IsHeadstageAvailableEvent [add, remove, raise]

 Event fires when 'true' if the headstage is connected for the headstage to query has changed

Additional Inherited Members

11.105.1 Detailed Description

CSCUFunctionNet is the class to control the SCU device

11.105.2 Constructor & Destructor Documentation

```
11.105.2.1 CSCUFunctionNet() [1/2] CSCUFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pSCUFunctionPointerContainer)
```

Initializes a new instance of the CSCUFunctionNet class.

```
11.105.2.2 CSCUFunctionNet() [2/2] CSCUFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.105.2.3 ~CSCUFunctionNet() virtual ~CSCUFunctionNet ( ) [virtual]
```

```
11.105.2.4 "!CSCUFunctionNet() !CSCUFunctionNet ()
```

11.105.3 Member Function Documentation

Sets automatic source channel selection

```
11.105.3.2 EnableAnalogOut() void EnableAnalogOut (
```

bool enable)

Enables AnalogOut globally

Parameters

```
enable Enable
```

11.105.3.3 GetAnalogOutADCRange() uint32_t GetAnalogOutADCRange ()

Gets the analog out ADC range

Returns

Range

Get the connected source channel number for a certain output channel

Parameters

out_channel	Output channel number

Returns

Source channel number

11.105.3.5 GetAnalogOutDACRange() AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ()

Gets the analog out DAC range

Returns

Range

$\textbf{11.105.3.6} \quad \textbf{GetAvailableHeadstages()} \quad \texttt{uint32_t GetAvailableHeadstages ()}$

Gets a bitmap of available headstages.

Returns

The bitmap of available headstages.

Gets multiple filter properties

Parameters

GroupID	The group ID
filterConfigurations_Length	The maximal length of filterConfigurations.

Returns

array of filter properties

Gets the filter property

Parameters

GroupID	The group ID
FilterNumber	The filter number

Returns

The filter property

```
11.105.3.9 GetHeadstageAdcBits() uint32_t GetHeadstageAdcBits ( uint32_t Headstage )
```

Gets the Number of ADC bits for a given headstage.

Headstage The headstage to query.

Returns

The number of bits the ADC has for the given headstage.

11.105.3.10 GetHeadstageAdcRangeInMicroVolt() uint32_t GetHeadstageAdcRangeInMicroVolt (uint32_t Headstage)

Gets the ADC Range in uV for a given headstage.

Parameters

Returns

The ADC Range in uV for the given headstage.

11.105.3.11 GetHeadstageDacBits() uint32_t GetHeadstageDacBits (uint32_t Headstage)

Gets the Number of DAC bits for a given headstage.

Parameters

Headstage	The headstage to query.
	, ,

Returns

The number of bits the DAC has for the given headstage.

Gets the DAC Current Range in uA for a given headstage.

Headstage	The headstage to query.
-----------	-------------------------

Returns

The DAC Current Range in uA for the given headstage.

11.105.3.13 GetHeadstageDacCurrentResolutionInNanoAmpere() uint32_t GetHeadstageDacCurrent↔

Gets the DAC Current Resolution in nA for a given headstage.

Parameters

Headstage	The headstage to query.
ouuotugo	in the measurage to query.

Returns

The DAC Current Resolution in nA for the given headstage.

Gets the DAC Voltage Range in mV for a given headstage.

Parameters

Headstage	The headstage to query.

Returns

The DAC Voltage Range in mV for the given headstage.

Gets the DAC Voltage Resolution in uV for a given headstage.

Headstage	The headstage to query.
-----------	-------------------------

Returns

The DAC Voltage Resolution in uV for the given headstage.

11.105.3.16 GetHeadstageGainInPermille() uint32_t GetHeadstageGainInPermille (uint32_t Headstage)

Gets the gain factor in permille for a given headstage.

Parameters

Returns

The gain factor in permille for the given headstage.

11.105.3.17 **GetHeadstageID()** SCU_HeadstageIdEnumNet GetHeadstageID (uint32_t Headstage)

Gets the headstage fpga ID.

Parameters

Headstane	The headstage to query.
Heaustage	The headstage to query.

Returns

The headstage fpga ID.

Gets the number of analog channels for a given headstage.

Returns

The number of analog channels the headstage has.

Gets the number of stimulation channels for a given headstage.

Parameters

Headstage The headstage to

Returns

The number of stimulation channels the headstage has.

11.105.3.20 GetHeadstageSamplerate() uint32_t GetHeadstageSamplerate (uint32_t Headstage)

Gets the Samplerate of a given headstage.

Parameters

Headstage	The headstage to query.

Returns

The samplerate in Hz for the given headstage.

```
11.105.3.21 GetHeadstageSerialNumber() String ^{\land} GetHeadstageSerialNumber ( uint32_t Headstage)
```

Gets the serial number of a given headstage.

Headstage	The headstage to query.
-----------	-------------------------

Returns

The serial number of the headstage.

11.105.3.22 GetMaxNumberOfHeadstages() uint32_t GetMaxNumberOfHeadstages ()

Gets the maximal number of headstages.

Returns

The maximal number of headstages.

11.105.3.23 GetMaxStimulusChannelsPerHeadstage() uint32_t GetMaxStimulusChannelsPerHeadstage

Gets the maximal number of stimulation channels a headstage can have.

Returns

The maximal number of stimulation channels a headstage can have.

Gets the mode for the reference electrode

Parameters

Headstage	The headstage number
-----------	----------------------

Returns

The mode

```
\textbf{11.105.3.25} \quad \textbf{GetReferenceElectrodeSwitchState()} \quad \texttt{ReferenceElectrodeSwitchPositionEnumNet} \quad \texttt{Get} \leftarrow \textbf{Continue} \quad \texttt{Continue} \quad \texttt
ReferenceElectrodeSwitchState (
                                                                                              uint32_t Headstage )
 Gets the position of the switch for the reference electrode
  Parameters
           Headstage
                                                                                               The headstage number
 Returns
                                   The switch position
11.105.3.26 HasAnalogOut() bool HasAnalogOut ( )
Has AnalogOut hardware
 Returns
                                   Enabled
11.105.3.27 HasGalvanicIsolation() bool HasGalvanicIsolation ( )
Has galvanic isolated hardware
 Returns
                                   Enabled
11.105.3.28 HasHSPowerSwitch() bool HasHSPowerSwitch ( )
Has SCU HS power switch
 Returns
                                   Has Switch
```

```
11.105.3.29 IsAnalogOutEnabled() bool IsAnalogOutEnabled ( )
```

Is AnalogOut enabled

Returns

Enabled

11.105.3.30 IsAutomaticAnalogOut() bool IsAutomaticAnalogOut ()

Is Automatic source channel selection selected

Returns

Automatic

```
11.105.3.31 IsHeadstageAvailable() bool IsHeadstageAvailable ( uint32_t Headstage )
```

Checks whether the given headstage is available.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

'true' if the headstage is connected.

```
11.105.3.32 IsHSPowered() bool IsHSPowered ( uint32_t Headstage )
```

Is the HS powered

Parameters

Headstage	The headstage to query.

Returns

'true' if the headstage is powered.

```
11.105.3.33 | IsInDacqLegacyMode() | bool | IsInDacqLegacyMode ( )
```

Is the SCU in legacy mode

Returns

Is Enabled

```
11.105.3.34 OnGetAvailableHeadstages() delegate void OnGetAvailableHeadstages ( uint32_t AvailableHeadstages )
```

```
11.105.3.35 OnlsHeadstageAvailable() delegate void OnIsHeadstageAvailable ( uint32_t Headstage, bool available )
```

Power the HS

Parameters

Headstage	The headstage to query.
power	'true' if the headstage is powered.

```
11.105.3.37 SetAnalogOutADCRange() void SetAnalogOutADCRange ( uint32_t range )
```

Sets the analog out ADC range

Parameters

range Range

11.105 CSCUFunctionNet Class Reference 434 Set the source channel number for a certain output channel

out_channel	Output channel number
source_channel	Source channel number

11.105.3.39 SetAnalogOutDACRange() void SetAnalogOutDACRange (

AnalogOut_DAC_Range_EnumNet range)

Sets the analog out DAC range

Parameters

range i

$\textbf{11.105.3.40} \quad \textbf{SetDacqLegacyMode()} \quad \texttt{void SetDacqLegacyMode ()}$

bool enable)

Enable the SCU legacy mode

Parameters

enable Enable

$\textbf{11.105.3.41} \quad \textbf{SetReferenceElectrodeMode()} \quad \texttt{void SetReferenceElectrodeMode} \quad \textbf{(}$

uint32_t Headstage,
ReferenceElectrodeModeEnumNet NewValue)

Sets the mode for the reference electrode

Parameters

Headstage	The headstage number
NewValue	The mode

$\textbf{11.105.3.42} \quad \textbf{SetReferenceElectrodeSwitchState()} \quad \texttt{void SetReferenceElectrodeSwitchState} \quad \textbf{(}$

uint32_t Headstage,

 ${\tt ReferenceElectrodeSwitchPositionEnumNet}\ \textit{NewSwitchPos}\)$

Sets the position of the switch for the reference electrode

Headstage	The headstage number
NewSwitchPos	The switch position

11.105.4 Event Documentation

11.105.4.1 GetAvailableHeadstagesEvent OnGetAvailableHeadstages^ GetAvailableHeadstagesEvent [add], [remove], [raise]


Event fires when the bitmap of available headstages has changed

11.105.4.2 IsHeadstageAvailableEvent OnIsHeadstageAvailable^ IsHeadstageAvailableEvent [add], [remove], [raise]

Event fires when 'true' if the headstage is connected for the headstage to query has changed

11.106 CSerialPortNet Class Reference

Inheritance diagram for CSerialPortNet:

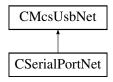

Public Member Functions

- CSerialPortNet (void)
- void Send (array< byte >^ buffer)
- void Send (String[^] command)
- array< byte > ^ Receive (void)
- array< byte > ^ Receive (int length)
- String ^ ReceiveString (void)
- String ^ ReceiveString (int length)
- int GetBytesAvailable (void)

Additional Inherited Members


11.106.1 Constructor & Destructor Documentation

```
11.106.1.1 CSerialPortNet() CSerialPortNet (
             void )
11.106.2 Member Function Documentation
11.106.2.1 GetBytesAvailable() int GetBytesAvailable (
             void )
11.106.2.2 Receive() [1/2] array<byte> ^ Receive (
             int length )
11.106.2.3 Receive() [2/2] array<br/>byte> ^{\wedge} Receive (
             void )
11.106.2.4 ReceiveString() [1/2] String ^ ReceiveString (
             int length )
11.106.2.5 ReceiveString() [2/2] String ^{\wedge} ReceiveString (
             void )
11.106.2.6 Send() [1/2] void Send (
             array< byte >^{\land} buffer )
11.106.2.7 Send() [2/2] void Send (
             String^{\wedge} command)
```

11.107 CStg200xBasicNet Class Reference

Base class for the Stg200x.

Inheritance diagram for CStg200xBasicNet:

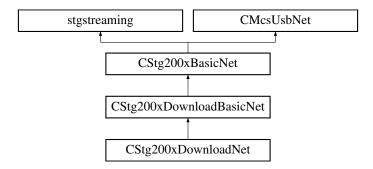

Public Member Functions

virtual ∼CStg200xBasicNet ()

The destructor.

void SetOutputRate (uint32_t rate)

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

• uint32 t GetOutputRate ()

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

void SendStart (uint32_t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

void SendStop (uint32_t triggermap)

Stop some or all triggers of the STG.

void SendStop (uint32_t triggermap, int options)

Stop some or all triggers of the STG.

• void GetStgVersionInfo ([Out]String^% SwVersion, [Out]String^% HwVersion)

Queries software and hardware version.

void GetAnalogRanges (int channel, [Out]int% URange, [Out]int% IRange)

Gets the range of the analog outputs.

• void GetAnalogResolution (int channel, [Out]int% URes, [Out]int% IRes)

Gets the resolution of the analog outputs.

virtual int32_t GetDACResolution ()

Gets number of bits of the DAC resolution.

virtual int32_t GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

virtual int32 t GetVoltageResolutionInMicroVolt (uint32 t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

• virtual int32_t GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

virtual int32_t GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

void GetStgProgramInfo ([Out]bool% IsProgrammed, [Out]System::Runtime::InteropServices::ComTypes::←
FILETIME% timestamp, [Out]String^% filename, [Out]Guid% guid)

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

void GetStgProgramInfo ([Out]bool% IsProgrammed, [Out]DateTime% timestamp, [Out]String^% filename, [Out]Guid% guid)

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

void SetStgProgramInfo (DateTime timestamp, String[^] filename, Guid guid)

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

• uint32 t GetAvailableMemory ()

Gets the amount of memory available in the currently selected segment of the STG.

uint32 t GetTotalMemory ()

Gets the total amount of memory available on the STG (all segments).

virtual uint32 t GetNumberOfAnalogChannels ()

Gets the Number of available analog channels of the device.

virtual uint32 t GetNumberOfSyncoutChannels ()

Gets the Number of available syncout channels of the device.

• virtual uint32_t GetNumberOfTriggerInputs ()

Gets the Number of trigger inputs of the device.

virtual uint32 t GetNumberOfHWDACPaths ()

Gets the Number of HW Stimulation DACs of the device.

virtual uint32 t GetNumberOfStimulationSourcesPerElectrode ()

Gets the number of stimulation sources (DACs) per electrode.

virtual void SetVoltageMode (unsigned int channel)

Sets a channel to voltage mode (STG3008-FA and STG400x only).

· virtual void SetCurrentMode (unsigned int channel)

Sets a channel to current mode (STG3008-FA and STG400x only).

virtual void SetVoltageMode ()

Sets all channels to voltage mode (STG3008-FA and STG400x only).

virtual void SetCurrentMode ()

Sets all channels to current mode (STG3008-FA and STG400x only).

virtual void SetMeasurementMode (unsigned int channel)

Sets a channel to measurement mode (STG3008-FA).

- · virtual void SetFAAmplification (unsigned int amplification)
- · virtual uint32 t GetFAAmplification ()
- virtual void SetAutocalibrationDisabled (unsigned int channel, bool disable)

Sets the autocalibration configuration.

virtual bool GetAutocalibrationDisabled (unsigned int channel)

Gets the autocalibration configuration.

virtual void SetElectrodeMode (uint32 t electrode, array< ElectrodeModeEnumNet >^ mode)

Puts an electrode in either automatic or manual mode.

• virtual void SetElectrodeMode (uint32_t electrode, ElectrodeModeEnumNet mode)

Puts an electrode in either automatic or manual mode.

virtual void SetElectrodeMode (uint32_t Scu_HS, uint32_t electrode, array< ElectrodeModeEnumNet >^ mode)

Puts an electrode in either automatic or manual mode.

• virtual void SetElectrodeMode (uint32_t Scu_HS, uint32_t electrode, ElectrodeModeEnumNet mode)

Puts an electrode in either automatic or manual mode.

• virtual uint32_t GetElectrodeMode (uint32_t electrode)

Gets the mode an electrode is in.

• virtual uint32_t GetElectrodeMode (uint32_t Scu_HS, uint32_t electrode)

Gets the mode an electrode is in.

Defines the DAC to use for an electrode.

virtual void SetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex, ElectrodeDacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

virtual void SetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, Electrode
 —
 DacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

virtual void SetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, array
 ElectrodeDacMuxEnumNet >^ dacMux)

Defines the DAC to use for an electrode.

virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex)

Gets the DAC which is used for an electrode.

virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_
 t listmodeIndex)

Gets the DAC which is used for an electrode.

virtual void SetElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an electrode.

virtual void SetElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

• virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, array< bool
 <p>>^ enable)

Enables or disables the stimulation switch for an electrode.

virtual bool GetElectrodeEnable (uint32_t electrode, uint32_t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

virtual bool GetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

virtual void SetExternalElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an external electrode.

virtual void SetExternalElectrodeEnable (uint32_t electrode, uint32_t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an external electrode.

virtual bool GetExternalElectrodeEnable (uint32_t electrode, uint32_t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

virtual void SetBlankingEnable (uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual void SetBlankingEnable (uint32_t electrode, array< bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

• virtual void SetBlankingEnable (uint32_t Scu_HS, uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual void SetBlankingEnable (uint32_t Scu_HS, uint32_t electrode, array < bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual bool GetBlankingEnable (uint32 t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

virtual bool GetBlankingEnable (uint32 t Scu HS, uint32 t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32_t electrode, bool enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32_t electrode, array< bool >^ enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

• virtual void SetEnableAmplifierProtectionSwitch (uint32 t Scu HS, uint32 t electrode, bool enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual bool GetEnableAmplifierProtectionSwitch (uint32_t electrode)

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual bool GetEnableAmplifierProtectionSwitch (uint32 t Scu HS, uint32 t electrode)

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

- virtual uint32 t GetNumberOfStimulationElectrodes ()
- template<typename digitalsourceenum >
 virtual void SetTriggerSource (unsigned int triggernum, DigitalSource< digitalsourceenum >^ triggersource,
 int bitnum offset)
- virtual void SetTriggerSource (unsigned int triggernum, TriggerSourceEnumNet triggersource, int bitnum_

 offset)
- virtual void SetTriggerSource (unsigned int triggernum, TriggerSourceEnumNet triggersource)
- virtual TriggerSourceEnumNet GetTriggerSource (unsigned int triggernum)
- virtual void SetListmodeIndexRange (unsigned int sideband, unsigned int startIndex, unsigned int endIndex, unsigned int mode)
- virtual void GetListmodeIndexRange (unsigned int sideband, unsigned int &startIndex, unsigned int &end← Index, unsigned int &mode)
- virtual void SetListmodeTriggerSource (unsigned int sideband, TriggerSourceEnumNet triggersource, int bitnumOffset)
- virtual void SetListmodeTriggerSource (unsigned int sideband, TriggerSourceEnumNet triggersource)
- virtual TriggerSourceEnumNet GetListmodeTriggerSource (unsigned int sideband)
- virtual void ListModeSendStart (unsigned int sidebandMask)
- virtual void ListModeSendStop (unsigned int sidebandMask)
- virtual void SetHeadstage (unsigned int headstage)
- virtual uint32_t GetHeadstage ()
- virtual void SetDacAmplificationFactor (uint32 t DacNumber, double Factor)

Set the amplification factor for a DAC.

virtual double GetDacAmplificationFactor (uint32 t DacNumber)

Get the amplification factor for a DAC.

• virtual void SetDigoutMode (Stg200xDigoutModeEnumNet digoutMode)

Sets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

virtual Stg200xDigoutModeEnumNet GetDigoutMode ()

Gets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

virtual void SetDigoutValue (uint32_t digoutValue)

Sets the Value on the digital output port when in manual mode.

virtual uint32 t GetDigoutValue ()

Gets the Value on the digital output port.

virtual uint32_t GetDiginValue ()

Gets the Value on the digital input port.

virtual void SetSyncoutMap (uint32_t channel, uint32_t syncoutMap)

Sets the mapping between external syncout outputs and internal syncout channels.

virtual uint32_t GetSyncoutMap (uint32_t channel)

Gets the mapping between external syncout outputs and internal syncout channels.

Additional Inherited Members

11.107.1 Detailed Description

Base class for the Stg200x.

From this class all STG related classes are derived: Mcs.Usb.CStg200xDownloadBasicNet Mcs.Usb.CStg200xDownloadNet for Download Mode and Mcs.Usb.CStg200xStreamingNet for Streaming Mode.

CStg200xBasicNet is the base class to control MCS STG device.

11.107.2 Constructor & Destructor Documentation

```
11.107.2.1 ~CStg200xBasicNet() virtual ~CStg200xBasicNet () [virtual]
```

The destructor.

11.107.3 Member Function Documentation

Gets the range of the analog outputs.

channel	The channel which is queried.
URange	The Voltage range in mV.
IRange	The Current range in uA.

```
11.107.3.2 GetAnalogResolution() void GetAnalogResolution ( int channel,
```

```
[Out] int% URes,
[Out] int% IRes )
```

Gets the resolution of the analog outputs.

Parameters

channel	The channel which is queried.
URes	The Voltage resolution in mV.
IRes	The Current resolution in uA.

11.107.3.3 **GetAutocalibrationDisabled()** virtual bool GetAutocalibrationDisabled (unsigned int *channel*) [virtual]

Gets the autocalibration configuration.

Parameters

channel The channel number.

Returns

true if autocalibration is disabled.

11.107.3.4 GetAvailableMemory() uint32_t GetAvailableMemory ()

Gets the amount of memory available in the currently selected segment of the STG.

Returns

The memory available in the currently selected segment in bytes.

```
11.107.3.5 GetBlankingEnable() [1/2] virtual bool GetBlankingEnable ( uint32_t electrode ) [virtual]
```

Gets whether an electrode should be blanked while stimulation is in progress.

electrode	The electrode number.
electroae	The electrode number.

Returns

true if blanking is enabled while stimulation is in progress.

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

Parameters

	electrode	The electrode number.
--	-----------	-----------------------

Returns

true if blanking is enabled while stimulation is in progress.

11.107.3.7 **GetCurrentRangeInNanoAmp()** virtual int32_t GetCurrentRangeInNanoAmp (uint32_t channel) [virtual]

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel Channel which is queried.	Channel which is queried	channel
-----------------------------------	--------------------------	---------

Returns

The Current Range of the specified channel in Nanoamps.

```
11.107.3.8 GetCurrentResolutionInNanoAmp() virtual int32_t GetCurrentResolutionInNanoAmp ( uint32_t channel) [virtual]
```

Gets the Current Resolution of the specified channel in Nanoamps.

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

```
11.107.3.9 GetDacAmplificationFactor() virtual double GetDacAmplificationFactor ( uint32_t DacNumber ) [virtual]
```

Get the amplification factor for a DAC.

Parameters

DacNumber The number of the DAC

Returns

the amplification factor for the DAC queried, range is from -1.99999 to +1.99999.

11.107.3.10 GetDACResolution() virtual int32_t GetDACResolution () [virtual]

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.107.3.11 GetDiginValue() virtual uint32_t GetDiginValue () [virtual]

Gets the Value on the digital input port.

Returns

The current value on the digital inputs.

11.107.3.12 GetDigoutMode() virtual Stg200xDigoutModeEnumNet GetDigoutMode () [virtual]

Gets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

Returns

The current operation mode.

11.107.3.13 GetDigoutValue() virtual uint32_t GetDigoutValue () [virtual]

Gets the Value on the digital output port.

Returns

The current value on the digital outputs.

Gets the DAC which is used for an electrode.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

11.107.3.15 GetElectrodeDacMux() [2/2] virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex) [virtual]

Gets the DAC which is used for an electrode.

Scu HS	The SCU headstage number.
0000	ine eee meaderage manneen

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

Scu HS	The SCU headstage number.
	· ·

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

```
11.107.3.18 GetElectrodeMode() [1/2] virtual uint32_t GetElectrodeMode ( uint32_t electrode ) [virtual]
```

Gets the mode an electrode is in.

Parameters

electrode	The electrode number.
Cicotioac	i i i i i i i i i i i i i i i i i i i

Returns

0 for automatic and 3 for manual mode.

```
11.107.3.19 GetElectrodeMode() [2/2] virtual uint32_t GetElectrodeMode ( uint32_t Scu_HS, uint32_t electrode ) [virtual]
```

Gets the mode an electrode is in.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
-----------	-----------------------

Returns

0 for automatic and 3 for manual mode.

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

The electrode numb	electrode
--------------------	-----------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

Gets weather an electrode is enabled or disabled for stimulation.

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

```
11.107.3.23 GetFAAmplification() virtual uint32_t GetFAAmplification ( ) [virtual]
```

```
11.107.3.24 GetHeadstage() virtual uint32_t GetHeadstage ( ) [virtual]
```

```
\textbf{11.107.3.25} \quad \textbf{GetListmodeIndexRange()} \quad \texttt{virtual void GetListmodeIndexRange ()}
```

```
unsigned int sideband,
unsigned int & startIndex,
unsigned int & endIndex,
unsigned int & mode ) [virtual]
```

```
11.107.3.26 GetListmodeTriggerSource() virtual TriggerSourceEnumNet GetListmodeTriggerSource ( unsigned int sideband ) [virtual]
```

11.107.3.27 GetNumberOfAnalogChannels() virtual uint32_t GetNumberOfAnalogChannels () [virtual]

Gets the Number of available analog channels of the device.

Returns

The number of analog channels.

11.107.3.28 GetNumberOfHWDACPaths() virtual uint32_t GetNumberOfHWDACPaths () [virtual]

Gets the Number of HW Stimulation DACs of the device.

Returns

The number of independent HW Stimulation outputs.

```
11.107.3.29 GetNumberOfStimulationElectrodes() virtual uint32_t GetNumberOfStimulationElectrodes ( ) [virtual]
```

```
11.107.3.30 GetNumberOfStimulationSourcesPerElectrode() virtual uint32_t GetNumberOfStimulation← SourcesPerElectrode ( ) [virtual]
```

Gets the number of stimulation sources (DACs) per electrode.

Returns

The number of stimulation sources (DACs) per electrode.

```
11.107.3.31 GetNumberOfSyncoutChannels() virtual uint32_t GetNumberOfSyncoutChannels ( ) [virtual]
```

Gets the Number of available syncout channels of the device.

Returns

The number of analog channels.

11.107.3.32 GetNumberOfTriggerInputs() virtual uint32_t GetNumberOfTriggerInputs () [virtual] Gets the Number of trigger inputs of the device.

Returns

The number of trigger inputs.

```
11.107.3.33 GetOutputRate() uint32_t GetOutputRate ( )
```

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Returns

Returns the current output rate in Hz.

11.107.3.34 GetStgProgramInfo() [1/2] void GetStgProgramInfo (

```
[Out] bool% IsProgrammed,
[Out] DateTime% timestamp,
[Out] String^% filename,
[Out] Guid% guid )
```

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

IsProgrammed	Flag wether download information is valid.
timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.
guid	A GUID.

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

Parameters

IsProgrammed	Flag wether download information is valid.
timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.

Queries software and hardware version.

Parameters

SwVersion	The current Software Version of the STG.
HwVersion	The Hardware Revision of the STG.

Gets the mapping between external syncout outputs and internal syncout channels.

channel	The external syncout output channel number (zero based).

Returns

The bitmap of internal syncout channels mapped to channel.

11.107.3.38 GetTotalMemory() uint32_t GetTotalMemory ()

Gets the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

```
11.107.3.39 GetTriggerSource() virtual TriggerSourceEnumNet GetTriggerSource ( unsigned int triggernum ) [virtual]
```

```
11.107.3.40 GetVoltageRangeInMicroVolt() virtual int32_t GetVoltageRangeInMicroVolt ( uint32_t channel) [virtual]
```

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

channel	Channel which is queried.

Returns

The Voltage Range of the specified channel in Microvolts.

```
11.107.3.41 GetVoltageResolutionInMicroVolt() virtual int32_t GetVoltageResolutionInMicroVolt ( uint32_t channel ) [virtual]
```

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

-		
ſ	channel	Channel which is queried.

Returns

The Voltage Resolution of the specified channel in Microvolts.

```
11.107.3.42 ListModeSendStart() virtual void ListModeSendStart ( unsigned int sidebandMask ) [virtual]
```

```
11.107.3.43 ListModeSendStop() virtual void ListModeSendStop ( unsigned int sidebandMask ) [virtual]
```

```
11.107.3.44 SendStart() void SendStart (
uint32_t triggermap )
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

triggermap	A bitmap of triggers which will be started.
------------	---

```
11.107.3.45 SendStop() [1/2] void SendStop ( uint32_t triggermap )
```

Stop some or all triggers of the STG.

Parameters

```
11.107.3.46 SendStop() [2/2] void SendStop ( uint32_t triggermap, int options )
```

Stop some or all triggers of the STG.

triggermap	A bitmap of triggers which will be stopped.	
options	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout assossiated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active	

11.107.3.47 SetAutocalibrationDisabled() virtual void SetAutocalibrationDisabled (unsigned int channel, bool disable) [virtual]

Sets the autocalibration configuration.

Parameters

channel	The channel number.
disable	true if autocalibration is to be disabled.

```
11.107.3.48 SetBlankingEnable() [1/4] virtual void SetBlankingEnable ( uint32_t electrode, array< bool >^ enable ) [virtual]
```

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether an electrode should be blanked while stimulation is in progress.

Scu_HS	The SCU headstage number.
--------	---------------------------

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

```
11.107.3.52 SetCurrentMode() [1/2] virtual void SetCurrentMode ( ) [virtual]
```

Sets all channels to current mode (STG3008-FA and STG400x only).

```
11.107.3.53 SetCurrentMode() [2/2] virtual void SetCurrentMode ( unsigned int channel ) [virtual]
```

Sets a channel to current mode (STG3008-FA and STG400x only).

channel	The channel to change.
---------	------------------------

Set the amplification factor for a DAC.

Parameters

DacNumber	The number of the DAC.
Factor	the amplification factor for that DAC, range is from -1.99999 to +1.99999.

11.107.3.55 SetDigoutMode() virtual void SetDigoutMode (Stg200xDigoutModeEnumNet digoutMode) [virtual]

Sets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

Parameters

	digoutMode	The new operation mode.
--	------------	-------------------------

```
11.107.3.56 SetDigoutValue() virtual void SetDigoutValue ( uint32_t digoutValue ) [virtual]
```

Sets the Value on the digital output port when in manual mode.

Parameters

digoutValue	The new value on the digital outputs.
-------------	---------------------------------------

Defines the DAC to use for an electrode.

electrode The elect	trode number.
---------------------	---------------

listmodeIndex	The index for listmode.	
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2	
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).	

Defines the DAC to use for an electrode.

Parameters

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use
	ElectrodeDacMuxEnumNet.Ground (0).

Defines the DAC to use for an electrode.

Scu_HS	The SCU headstage number.

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use
	ElectrodeDacMuxEnumNet.Ground (0).

$\textbf{11.107.3.60} \quad \textbf{SetElectrodeDacMux() [4/4]} \quad \text{virtual void SetElectrodeDacMux ()}$

```
uint32_t Scu_HS,
uint32_t electrode,
uint32_t listmodeIndex,
ElectrodeDacMuxEnumNet dacMux ) [virtual]
```

Defines the DAC to use for an electrode.

Parameters

Scu_HS The SCU headstage number.

Parameters

electrode	The electrode number.
electione	i i ne electrode number.

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use
	ElectrodeDacMuxEnumNet.Ground (0).

11.107.3.61 SetElectrodeEnable() [1/4] virtual void SetElectrodeEnable (uint32_t electrode,

```
uint32_t listmodeIndex, array < bool >^{ } enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an electrode.

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an electrode.

Scu_HS The SCU headstage number.	Scu_HS	The SCU headstage number.
------------------------------------	--------	---------------------------

Parameters

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

$\textbf{11.107.3.64} \quad \textbf{SetElectrodeEnable() [4/4]} \quad \text{virtual void SetElectrodeEnable (}$

```
uint32_t Scu_HS,
uint32_t electrode,
uint32_t listmodeIndex,
bool enable ) [virtual]
```

Enables or disables the stimulation switch for an electrode.

Parameters

Scu HS The SCU neadstage number.	Scu HS	The SCU headstage number.
------------------------------------	--------	---------------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

```
11.107.3.65 SetElectrodeMode() [1/4] virtual void SetElectrodeMode ( uint32_t electrode, array< ElectrodeModeEnumNet >^ mode ) [virtual]
```

Puts an electrode in either automatic or manual mode.

Parameters

electrode	The electrode number.
-----------	-----------------------

Returns

0 for automatic and 3 for manual mode.

Puts an electrode in either automatic or manual mode.

Parameters

electrode	The electrode number.
electione	I THE ELECTIONS HUTTIDET.

Parameters

mode	0 for automatic and 3 for manual mode.
IIIOUE	, o loi automatic and 3 loi manuai mode.

Puts an electrode in either automatic or manual mode.

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
-----------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.107.3.68 SetElectrodeMode() [4/4] virtual void SetElectrodeMode (

```
uint32_t Scu_HS,
uint32_t electrode,
ElectrodeModeEnumNet mode ) [virtual]
```

Puts an electrode in either automatic or manual mode.

Parameters

Parameters

Parameters

mode 0 for automatic and 3 for manual mode.

11.107.3.69 SetEnableAmplifierProtectionSwitch() [1/4] virtual void SetEnableAmplifierProtection←

```
uint32_t electrode,
array< bool >^{\wedge} enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

electrod	The electrode number.	1
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.	

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Scu_HS	The SCU headstage number.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.107.3.72 SetEnableAmplifierProtectionSwitch() [4/4] virtual void SetEnableAmplifierProtection← Switch (

```
uint32_t Scu_HS,
uint32_t electrode,
bool enable ) [virtual]
```

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Enables or disables the stimulation switch for an external electrode.

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an external electrode.

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

```
11.107.3.75 SetFAAmplification() virtual void SetFAAmplification (
             unsigned int amplification ) [virtual]
11.107.3.76 SetHeadstage() virtual void SetHeadstage (
             unsigned int headstage ) [virtual]
11.107.3.77 SetListmodeIndexRange() virtual void SetListmodeIndexRange (
             unsigned int sideband,
             unsigned int startIndex,
             unsigned int endIndex,
             unsigned int mode ) [virtual]
11.107.3.78 SetListmodeTriggerSource() [1/2] virtual void SetListmodeTriggerSource (
             unsigned int sideband,
             {\tt TriggerSourceEnumNet}\ triggersource\ )\quad [{\tt virtual}]
11.107.3.79 SetListmodeTriggerSource() [2/2] virtual void SetListmodeTriggerSource (
             unsigned int sideband,
             TriggerSourceEnumNet triggersource,
             int bitnumOffset ) [virtual]
11.107.3.80 SetMeasurementMode() virtual void SetMeasurementMode (
```

unsigned int channel) [virtual]

Sets a channel to measurement mode (STG3008-FA).

Generated by Doxygen

channel	The channel to change.
---------	------------------------

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Parameters

rate	The new output rate in Hz.
------	----------------------------

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

Parameters

timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.

Sets the mapping between external syncout outputs and internal syncout channels.

channel	The external syncout output channel number (zero based).
syncoutMap	A bitmap of internal syncout channels to map to channel.

```
11.107.3.84 SetTriggerSource() [1/3] virtual void SetTriggerSource ( unsigned int triggernum,
```

```
\label{eq:digitalsource} \begin{tabular}{ll} DigitalSource < digitalsource enum >^ $$ triggersource, \\ int $bitnum\_offset $) $$ [virtual] $$ \end{tabular}
```

```
11.107.3.85 SetTriggerSource() [2/3] virtual void SetTriggerSource ( unsigned int triggernum,
```

TriggerSourceEnumNet triggersource) [virtual]

11.107.3.86 SetTriggerSource() [3/3] virtual void SetTriggerSource (

```
unsigned int triggernum,
TriggerSourceEnumNet triggersource,
int bitnum_offset ) [virtual]
```

11.107.3.87 SetVoltageMode() [1/2] virtual void SetVoltageMode () [virtual]

Sets all channels to voltage mode (STG3008-FA and STG400x only).

```
11.107.3.88 SetVoltageMode() [2/2] virtual void SetVoltageMode (
unsigned int channel) [virtual]
```

Sets a channel to voltage mode (STG3008-FA and STG400x only).


Parameters

channel	The channel to change.

11.108 CStg200xDownloadBasicNet Class Reference

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

Inheritance diagram for CStg200xDownloadBasicNet:

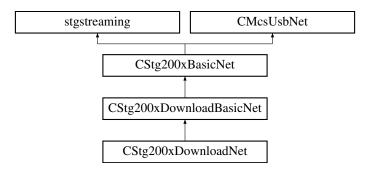

Public Member Functions

virtual void SetupTrigger (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

- virtual void SetupTriggerSingle (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)

 Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void GetTrigger ([Out] array< uint32_t >^% channelmap, [Out] array< uint32_t >^% syncoutmap, [Out] array< uint32_t >^% repeat)

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

void GetSweepCount ([Out] array< uint32 t >^% sweeps, [Out] array< uint32 t >^% triggers)

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in CStg200xDownloadBasicNet::SetupTrigger.
- void ForceStatusEvent ()

Force a status event.

void ResetStatus (uint32 t triggermap)

Reset the status flag.

uint32 t GetMemoryUsageDAC (uint32 t Channel)

Queries the memory usage of the current segment and selected analog DAC channel.

uint32_t GetMemoryUsageSyncout (uint32_t Channel)

Queries the memory usage of the current segment and selected syncout channel.

virtual void ClearSyncData (uint32_t channel)

Delete a SyncOut pattern for a channel from STG memory.

virtual void SendSyncData (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μs. The STG has a resolution of 20 μs. If your application can not handle 64 bit integers, use the STG200x SendSyncData32() call instead.

virtual void ClearChannelData (uint32_t channel)

Delete a stimulus pattern for a channel from STG memory

virtual void SendChannelData (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μs. The STG has a resolution of 20 μs.

virtual void EnableAutoReset ()

Enable AutoReset of the STG Status.

• virtual void DisableAutoReset ()

Disable AutoReset of the STG Status.

 virtual void SetupRetriggerMode (int8_t trigger, RetriggerActionEnumNet same_trigger, RetriggerAction← EnumNet other trigger)

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal
- virtual void SetupRetriggerMode (RetriggerActionEnumNet same_trigger, RetriggerActionEnumNet other_
 trigger)

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Properties

• CStimulusFunctionNet[^] Stimulus [get]

Additional Inherited Members

11.108.1 Detailed Description

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

11.108.2 Member Function Documentation

```
11.108.2.1 ClearChannelData() virtual void ClearChannelData ( uint32_t channel) [virtual]
```

Delete a stimulus pattern for a channel from STG memory

Parameters

channel Specifies the channel to clear.

```
11.108.2.2 ClearSyncData() virtual void ClearSyncData ( uint32_t channel ) [virtual]
```

Delete a SyncOut pattern for a channel from STG memory.

Parameters

channel Specifies the syncout channel to clear.

11.108.2.3 DisableAutoReset() virtual void DisableAutoReset () [virtual]

Disable AutoReset of the STG Status.

If autoreset is disabled, the STG status switches to FINISHED after the defined number of sweeps is finished. To switch back to the IDLE status, use CStg200xDownload::ResetStatus()

```
11.108.2.4 EnableAutoReset() virtual void EnableAutoReset ( ) [virtual]
```

Enable AutoReset of the STG Status.

This is the default on power up. If autoreset is enabled, the STG status switches to FINISHED only for one poll cycle after this, it switches to IDLE automatically.

```
11.108.2.5 ForceStatusEvent() void ForceStatusEvent ()
```

Force a status event.

Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

```
11.108.2.6 GetMemoryUsageDAC() uint32_t GetMemoryUsageDAC ( uint32_t Channel)
```

Queries the memory usage of the current segment and selected analog DAC channel.

The currently used memory is reported for the requested channel.

Parameters

Returns

Returns the usage in STG memory.

```
11.108.2.7 GetMemoryUsageSyncout() uint32_t GetMemoryUsageSyncout ( uint32_t Channel)
```

Queries the memory usage of the current segment and selected syncout channel.

The currently used memory is reported for the requested channel.

Channel	channel for the amount of interested usage.
---------	---

Returns

Returns the usage in STG memory.

```
11.108.2.8 GetSweepCount() void GetSweepCount ( [Out] \ array < \ uint 32\_t > ^{ } \% \ sweeps, [Out] \ array < \ uint 32\_t > ^{ } \% \ triggers )
```

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in CStg200xDownloadBasicNet::SetupTrigger.

Parameters

sweeps	on return contains the number of sweeps for each trigger.
triggers	on return contains the number of trigger events seen for each trigger.

```
11.108.2.9 GetTrigger() void GetTrigger (

[Out] array< uint32_t >^% channelmap,

[Out] array< uint32_t >^% syncoutmap,

[Out] array< uint32_t >^% repeat )
```

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

channelmap	For each trigger, a bitmap of channels that belong to this trigger.
	99-,

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

```
11.108.2.10 ResetStatus() void ResetStatus (
uint32_t triggermap )
```

Reset the status flag.

triggerm	nap bitma	o of trigger fo	or which to r	reset the status.
----------	-----------	-----------------	---------------	-------------------

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Parameters

channel	Specifies the channel to append the data to.	
pData	A list of datapoints.	
tData	A list of durations as int64_t. The time is given in units of μs.	

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the STG200x_SendSyncData32() call instead.

channel	Specifies the sync output channel to append the data to.		
pData	Data A list of datapoints.		
tData	A list of durations as int64_t. The time is given in units of μs.		

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- · restart this trigger
- · ignore the signal

Parameters

trigger	The trigger to change.
same_trigger	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
other_trigger	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected.Not used in Normal Mode.

```
11.108.2.14 SetupRetriggerMode() [2/2] virtual void SetupRetriggerMode (
RetriggerActionEnumNet same_trigger,
RetriggerActionEnumNet other_trigger) [virtual]
```

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- · restart this trigger
- · ignore the signal

Parameters

same_trigger	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
other_trigger	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

11.108.2.15 SetupTrigger() virtual void SetupTrigger (

```
uint32_t first_trigger,
array< uint32_t >^ channelmap,
array< uint32_t >^ syncoutmap,
array< uint32_t >^ repeat ) [virtual]
```

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

irst_trigger The number of the first trigger to cha	ige.
---	------

Parameters

channelmap	For each trigger, a bitmap of channels that belong to this trigger.
------------	---

Parameters

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

trigger	The trigger to change.

channelmap	A bitmap of channels that belong to this trigger.

syncoutmap	A bitmap of syncouts that belong to this trigger.
repeat	The number of times this trigger should be repeated.


11.108.3 Property Documentation

11.108.3.1 Stimulus CStimulusFunctionNet^ Stimulus [get]

11.109 CStg200xDownloadNet Class Reference

Main class for the STG download mode This class implements the STG download mode interface.

Inheritance diagram for CStg200xDownloadNet:

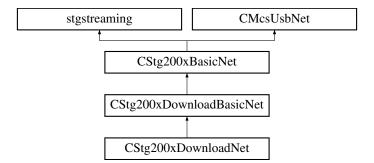

Public Member Functions

CStg200xDownloadNet ()

Use this constructor if you do not want to use the status callback.

CStg200xDownloadNet (OnStgPollStatus[^] pollStatus)

Use this constructor if you want to use the status callback.

- ∼CStg200xDownloadNet ()
- void PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

 void PrepareAndAppendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and append data to a given channel on the STG.

void ClearChannel_PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType, bool doClear)

Prepare and append data to a given channel on the STG.

void SendSegmentDefine (array< uint32_t >^ segment_list)

Defines the segment memory layout of the STG.

 void SendSegmentStart (uint32_t triggermap, uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)

Switchs segment and starts trigger.

- void SendSegmentSelect (uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)
 Switchs segment.
- void EnableMultiFileMode (uint32_t submode)

Enable the Multi-File mode of the STG.

void DisableMultiFileMode ()

Disable the Multi-File mode of the STG

- StgStatusNet ^ QueryTriggerstatus ()
- void SetOutputMap (array< uint32_t $>^{\land}$ ChannelLayout)
- int32 t GetModuleTemp (unsigned int channel)
- uint32_t GetModuleCurrent (unsigned int channel)

Events

- OnStgPollStatus Stg200xPollStatusEvent [add, remove, raise]
- OnMwPollStatus^ MwPollStatusEvent [add, remove, raise]

Additional Inherited Members

11.109.1 Detailed Description

Main class for the STG download mode This class implements the STG download mode interface.

11.109.2 Constructor & Destructor Documentation

```
11.109.2.1 CStg200xDownloadNet() [1/2] CStg200xDownloadNet ( )
```

Use this constructor if you do not want to use the status callback.

```
11.109.2.2 CStg200xDownloadNet() [2/2] CStg200xDownloadNet ( onStgPollStatus^{\land} pollStatus )
```

Use this constructor if you want to use the status callback.

```
11.109.2.3 ~CStg200xDownloadNet() ~CStg200xDownloadNet ( )
```

11.109.3 Member Function Documentation

Prepare and append data to a given channel on the STG.

bool doClear)

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel The channel number to send data	to.
---	-----

Parameters

amplitude A list of amplitudes in units of μV	and nA in voltage and current mode, respectively.
---	---

Parameters

duration	A list of durations in units of μs.
destType	specifies wheather the data is for syncout, current or voltage stimulation.

11.109.3.2 DisableMultiFileMode() void DisableMultiFileMode ()

Disable the Multi-File mode of the STG

Switch the STG back to normal mode. In this mode, trigger inputs are assigned to channels, not to segments.

```
11.109.3.3 EnableMultiFileMode() void EnableMultiFileMode ( uint32_t submode )
```

Enable the Multi-File mode of the STG.

In Multi-File mode, the trigger inputs switch between segments. To use this mode, define up to as many segments as trigger inputs are available and fill each segment with a stimulus pattern.

Now a trigger on trigger input 1 switches the STG to the first segment and starts all triggers in this segment. Likewise, a trigger on trigger input 2, 3 and 4 selects the respective segment and start all triggers in this segment So the Multi-File Mode can be used to predefine up to four different stimuli which can be selected without the need for a computer connection.

Parameters

submode

The submode. Submode 0 is regular Multi-File mode as described above, submode 1 is extended Multi-File mode, where the segment is selected based on the digital pattern on the digital inputs. In this mode, 256 different segments can be defined and used.

```
11.109.3.6 PrepareAndAppendData() void PrepareAndAppendData (
```

```
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType )
```

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

channel	The channel number to send data to.
---------	-------------------------------------

Parameters

amplitude A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
--

Parameters

duration	A list of durations in units of μs.
destType	specifies wheather the data is for syncout, current or voltage stimulation.

$\textbf{11.109.3.7} \quad \textbf{PrepareAndSendData()} \quad \texttt{void PrepareAndSendData} \ \ \textbf{(}$

```
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType )
```

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel	The channel number to send data to.
---------	-------------------------------------

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

duration	A list of durations in units of μs.	
destType specifies wheather the data is for syncout, current or voltage stimula		

11.109.3.8 QueryTriggerstatus() StgStatusNet ^ QueryTriggerstatus ()

```
11.109.3.9 SendSegmentDefine() void SendSegmentDefine ( array < uint32_t >^{\land} segment\_list )
```

Defines the segment memory layout of the STG.

On reset, the STG has one segment containing all available memory.

With this command, the STG memory can be devided into several segments. Each segment can be filled with stimulus data.

Parameters

	segment_list	The List of memory sizes (one per segment).
--	--------------	---

Switchs segment.

Parameters

segment	The number of the segment to select.

Switchs segment and starts trigger.

Parameters

Parameters

Parameters

segmentflags | A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.

```
11.109.3.12 SetOutputMap() void SetOutputMap (

array< uint32_t >^ ChannelLayout )
```


11.109.4 Event Documentation

```
11.109.4.1 MwPollStatusEvent OnMwPollStatus^ MwPollStatusEvent [add], [remove], [raise]
```

```
11.109.4.2 Stg200xPollStatusEvent OnStgPollStatus^ Stg200xPollStatusEvent [add], [remove], [raise]
```

11.110 CStimulusFunctionNet Class Reference

Inheritance diagram for CStimulusFunctionNet:

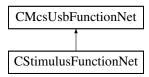

Classes

- class SidebandData
- class StimulusDeviceDataAndUnrolledData

Public Member Functions

- CStimulusFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] stimulusFunction
 —
 PointerContainer)
- CStimulusFunctionNet (CMcsUsbNet[^] mcsusb)
- void StartPoll ()

Starts the interrupt fetching thread and delivers events

void StopPoll ()

Stops the interrupt fetching thread and delivers events

void ForceStatusEvent ()

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

void SendStart (uint32_t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

void SendStop (uint32 t triggermap)

Stop some or all triggers of the STG.

void SendStop (uint32_t triggermap, int options)

Stop some or all triggers of the STG.

· void ClearChannelData (int channel)

Delete a Stimulus Pattern from STG memory

• void ClearSyncData (int channel)

Delete a Syncout Pattern from STG memory

 void PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

 void PrepareAndAppendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and append data to a given channel on the STG.

- void ClearChannel_PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType, bool doClear)
- StimulusDeviceDataAndUnrolledData ^ PrepareData (int channel, array< int32_t > ^ amplitude, array< uint64 t > ^ duration, STG DestinationEnumNet destType)
- void SendPreparedData (int channel, StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled, S
 —
 TG DestinationEnumNet destType)
- SidebandData ^ CreateSideband (array< int32_t >^ StimulusActive, array< int32_t >^ Syncout, array< uint64_t >^ Duration, uint32_t Bit0Time, uint32_t Bit3Time, uint32_t Bit4Time)

Creates the Sideband Channel for the MEA2100 device.

void ClearMultiplexedData ()

Clears the Stimulation Memory in the STG device.

void SendMultiplexedData (array< uint16_t >^ data)

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

• int GetMultiplexedDataChannelsInBlock ()

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData

• int GetDACResolution ()

Gets number of bits of the DAC resolution.

• int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

int GetVoltageResolutionInMicroVolt (uint32 t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

int GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

• int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

void SetupTrigger (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

void SetupTriggerSingle (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

uint32_t GetTotalMemory ()

Get the total amount of memory available on the STG (all segments).

uint32 t GetAvailableMemory ()

Get the amount of memory available in the currently selected segment of the STG.

int GetNumberOfAnalogChannels ()

Get the number of STG channels.

Events

OnStgPollStatus^ PollStatusEvent

Additional Inherited Members

11.110.1 Constructor & Destructor Documentation

```
11.110.1.1 CStimulusFunctionNet() [1/2] CStimulusFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ stimulusFunctionPointerContainer)
```

```
11.110.1.2 CStimulusFunctionNet() [2/2] CStimulusFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.110.2 Member Function Documentation

Delete a Stimulus Pattern from STG memory

Parameters

channel specifies the channel to clear.

11.110.2.3 ClearMultiplexedData() void ClearMultiplexedData ()

Clears the Stimulation Memory in the STG device.

Delete a Syncout Pattern from STG memory

Parameters

channel | specifies the channel to clear.

Creates the Sideband Channel for the MEA2100 device.

Each datapoint is represented by an signed 32bit integer value. A value 0 means that the stimulation is active during that time. A value 1 means that the stimulation is not active during that time.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Parameters

StimulusActive	A list of datapoints which define weather the Stimulus is active or idle at that time as int32.
----------------	---

Parameters

Duration	A list of durations as uint64. The time is given in units of μs.
Bit0Time	Time in µs for which Bit 0 (Blanking) is to be extended.

Parameters

Bit3Time	Time in µs for which Bit 3 (Stimulus Enable) is to be extended.
----------	---

Parameters

Bit4Time	Time in µs for which Bit 4 (Stimulus Selector) is to be extended.
----------	---

Returns

Error Status. 0 on success.

11.110.2.6 ForceStatusEvent() void ForceStatusEvent ()

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

11.110.2.7 GetAvailableMemory() uint32_t GetAvailableMemory ()

Get the amount of memory available in the currently selected segment of the STG.

Returns

The total memory available on the STG in bytes.

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.110.2.9 GetCurrentResolutionInNanoAmp() int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.110.2.10 GetDACResolution() int GetDACResolution ()

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.110.2.11 GetMultiplexedDataChannelsInBlock() int GetMultiplexedDataChannelsInBlock ()

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData

11.110.2.12 GetNumberOfAnalogChannels() int GetNumberOfAnalogChannels ()

Get the number of STG channels.

Returns

The number of STG channels.

11.110.2.13 GetTotalMemory() uint32_t GetTotalMemory ()

Get the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

11.110.2.14 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

11.110.2.15 **GetVoltageResolutionInMicroVolt()** int GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

channel Channel which is queried.	
-----------------------------------	--

Returns

The Voltage Resolution of the specified channel in Microvolts.

11.110.2.16 PrepareAndAppendData() void PrepareAndAppendData (

```
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType )
```

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

C	hannel	The channel number to send data to.
---	--------	-------------------------------------

Parameters

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

duration	A list of durations in units of μs.
destType	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

11.110.2.17 PrepareAndSendData() void PrepareAndSendData (

```
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType )
```

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel	The channel number to send data to.
---------	-------------------------------------

Parameters

Parameters

duration	A list of durations in units of μs.
destType	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

```
11.110.2.18 PrepareData() StimulusDeviceDataAndUnrolledData ^ PrepareData (
    int channel,
    array< int32_t >^ amplitude,
    array< uint64_t >^ duration,
    STG_DestinationEnumNet destType )
```

```
11.110.2.19 SendMultiplexedData() void SendMultiplexedData ( array < uint16_t >^{\wedge} data )
```

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

Parameters

```
data Array of data to be sent.
```

```
11.110.2.21 SendStart() void SendStart ( uint32_t triggermap )
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

triggermap A bitmap of triggers which will be started.

```
11.110.2.22 SendStop() [1/2] void SendStop ( uint32_t triggermap )
```

Stop some or all triggers of the STG.

Parameters

triggermap	A bitmap of triggers which will be stopped.

```
11.110.2.23 SendStop() [2/2] void SendStop (
    uint32_t triggermap,
    int options)
```

Stop some or all triggers of the STG.

Parameters

triggermap	A bitmap of triggers which will be stopped.	
options	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses the stop commands when a syncout assossiated with a given sync-out has bit 1 (0x02) set. Can be used e.g. to prevent a stop while a biphasic stimulation pulse is active	

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

first_trigger The number of the first trigger to	change.
---	---------

Parameters

-1	
cnanneiman	For each trigger, a bitmap of channels that belong to this trigger.
oaap	i or each ingger, a chinap or chamber in at coloning to time ingger

Parameters

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.	

Parameters

channelmap	A bitmap of channels that belong to this trigger.

Parameters

syncoutmap	A bitmap of syncouts that belong to this trigger.
repeat	The number of times this trigger should be repeated.

11.110.2.26 StartPoll() void StartPoll ()

Starts the interrupt fetching thread and delivers events

11.110.2.27 StopPoll() void StopPoll ()

Stops the interrupt fetching thread and delivers events


11.110.3 Event Documentation

11.110.3.1 PollStatusEvent OnStgPollStatus^ PollStatusEvent

11.111 CSw2to64DeviceNet Class Reference

The class to control the MCS-USB-Sw2to64 device.

Inheritance diagram for CSw2to64DeviceNet:

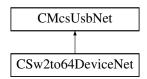

Public Member Functions

- CSw2to64DeviceNet ()
- ∼CSw2to64DeviceNet ()
- unsigned short GetNumber ()

Gets the number of channels that can be switched in this box.

array< unsigned char > ^ GetChannels ()

Gets the current switch positions as char array.

void SetChannels (array< unsigned char >^ pattern)

Sets the switch positions from a char array.

unsigned char GetChannel (unsigned short index)

Gets one current switch position.

void SetChannel (unsigned short index, unsigned char pattern)

Sets one switch position.

Additional Inherited Members

11.111.1 Detailed Description

The class to control the MCS-USB-Sw2to64 device.

This class controls the settings of the MCS-USB-Sw2to64. The box has two inputs for signals. Each of the 64 outputs can be connected to one of the input signals, could be held open or connected ground. Valid switch states are 0, 1, 2 or 3 for each of the settings.

11.111.2 Constructor & Destructor Documentation

```
11.111.2.1 CSw2to64DeviceNet() CSw2to64DeviceNet ()
```

```
11.111.2.2 ~CSw2to64DeviceNet() ~CSw2to64DeviceNet ()
```

11.111.3 Member Function Documentation

```
11.111.3.1 GetChannel() unsigned char GetChannel ( unsigned short index )
```

Gets one current switch position.

in index number of channel to read the switch po	osition from
--	--------------

Returns

switch position of desired channel

11.111.3.2 GetChannels() array<unsigned char> ^ GetChannels ()

Gets the current switch positions as char array.

Returns

array of char with the size of the number of channels, each char has the setting of a channel

11.111.3.3 GetNumber() unsigned short GetNumber ()

Gets the number of channels that can be switched in this box.

The box can have a different number of channels it can switch. Up to now usually 64 channels are returned

```
11.111.3.4 SetChannel() void SetChannel (
unsigned short index,
unsigned char pattern )
```

Sets one switch position.

Parameters

i	n	index	number of channel to write the switch position to
i	n	pattern	switch position of the channel

11.111.3.5 SetChannels() void SetChannels ($array < unsigned char >^{\wedge} pattern$)

Sets the switch positions from a char array.


Parameters

in	pattern	array of char with the size of the number of channels, each char has the setting of a channel

11.112 CTcxDeviceNet Class Reference

Class to control a Temperature Controller (TCX)

Inheritance diagram for CTcxDeviceNet:

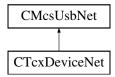

Public Member Functions

CTcxDeviceNet ()

Initializes a new instance of CTcxDeviceNet class.

- ∼CTcxDeviceNet ()
- unsigned int GetNumControlChannels ()

Gets the number of channels the device can control/regulate.

unsigned int GetNumMeasureChannels ()

Gets the number of channels the device can measure.

• int GetValue (unsigned int channel)

Gets the temperate of the specified channel in units of 0.1 ℃.

• int GetValueHires (unsigned int channel)

Gets the temperate of the specified channel in units of 0.01 ℃.

int GetHeaterTemp (unsigned int channel)

Gets the temperate of the specified heater in units of 0.1 °C.

int GetHeaterLimit (unsigned int device)

Gets the temperate limit of the specified heater in units of 0.1 °C.

double GetMaxHeaterPowerMultiwell ()

queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

void SetMaxHeaterPowerMultiwell (double MaxPowerWatt)

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

bool GetHasThermocouple ()

Gets weather the device supports a thermocouple.

- bool GetEnableHeaterLimit (unsigned int device)
- bool GetEnableThermocouple (unsigned int device)
- TcxSensorTypeEnumNet GetSensorType (unsigned int device)
- String \(^\) GetUnit (unsigned int channel)
- unsigned int GetBoardTemp ()

Gets the temperate of the mainboard in units of 0.1 °C.

- · unsigned int GetVolti (unsigned int channel)
- unsigned int GetNumDevices ()
- void SetSetpoint (unsigned int channel, int sp)

Sets the target temperate of specified channel in units of 0.1 °C.

- void SetDevice (unsigned int channel, int device)
- void SetOnOff (unsigned int channel, bool on)

Switches the specified channel on or off.

- · void SetCalibration (unsigned int channel, int calib)
- void SetP (unsigned int device, int p_coeff)

Sets the P-coefficient of the specified device.

void SetI (unsigned int device, int i coeff)

Sets the I-coefficient of the specified device.

void SetD (unsigned int device, int d_coeff)

Sets the D-coefficient of the specified device.

void SetMaxP (unsigned int device, int maxp)

Sets the maximum heater power of the specified device.

- void SetHeaterLimit (unsigned int device, int heater_limit)
- void SetEnableHeaterLimit (unsigned int device, bool enable)
- void SetEnableThermocouple (unsigned int device, bool enable)
- void SetSensorType (unsigned int device, TcxSensorTypeEnumNet type)
- void SetDevname (unsigned int device, String[^] Devicename)
- · int GetSetpoint (unsigned int channel)

Gets the target temperate of specified channel in units of 0.1 ℃.

- int GetDevice (unsigned int channel)
- int GetOnOff (unsigned int channel)

Gets if the specified channel is on or off.

- int GetCalibration (unsigned int channel)
- int GetP (unsigned int device)

Gets the P-coefficient of the specified device.

· int GetI (unsigned int device)

Gets the I-coefficient of the specified device.

int GetD (unsigned int device)

Gets the D-coefficient of the specified device.

• int GetMaxP (unsigned int device)

Gets the maximum heater power of the specified device.

- String \(^\) GetDevname (unsigned int device)
- TcxDeviceTypeEnumNet GetDeviceType ()
- int GetSetpointMin (unsigned int channel)
- int GetCalibrationMin (unsigned int channel)
- int GetPMin (unsigned int device)
- int GetIMin (unsigned int device)
- int GetDMin (unsigned int device)
- int GetMaxpMin (unsigned int device)
- · int GetSetpointMax (unsigned int channel)
- · int GetCalibrationMax (unsigned int channel)
- int GetPMax (unsigned int device)
- int GetIMax (unsigned int device)
- int GetDMax (unsigned int device)
- int GetMaxpMax (unsigned int device)
- int GetSetpointDecp (unsigned int channel)
- int GetCalibrationDecp (unsigned int channel)
- int GetPDecp (unsigned int device)
- int GetIDecp (unsigned int device)
- int GetDDecp (unsigned int device)
- int GetMaxpDecp (unsigned int device)
- · int GetResX (unsigned int channel)
- · int GetResS (unsigned int channel)
- int GetRes1 (unsigned int channel)
- int GetRes2 (unsigned int channel)
- int GetPwrSet (unsigned int channel)
- int GetPwrOut (unsigned int channel)
- int GetDuty (unsigned int channel)

Gets the duty cycle of the heating element.

• int GetUOut (unsigned int channel)

Gets the voltage on the heating element.

• int GetlOut (unsigned int channel)

Gets the current through the heating element.

int GetROut (unsigned int channel)

Gets the resistance of the heating element.

int GetPOut (unsigned int channel)

Gets the output power of the heating element.

- int GetCurrent (unsigned int channel)
- int GetThermocoupleTemp (unsigned int channel)
- int GetThermocoupleTempAbs (unsigned int channel)
- int GetThermocoupleReferenceTemp (unsigned int channel)
- unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)

Gets the proportional constant for the thermocouple.

void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)

Sets the proportional constant for the thermocouple.

- int GetThermocoupleCalibration (unsigned int channel)
- void CalibrateThermocouple (unsigned int channel)
- void SetDeviceType (TcxDeviceTypeEnumNet devicetype)
- void FactoryReset ()

Additional Inherited Members

11.112.1 Detailed Description

Class to control a Temperature Controller (TCX)

11.112.2 Constructor & Destructor Documentation

```
11.112.2.1 CTcxDeviceNet() CTcxDeviceNet ( )
```

Initializes a new instance of CTcxDeviceNet class.

```
11.112.2.2 ~CTcxDeviceNet() ~CTcxDeviceNet ()
```

11.112.3 Member Function Documentation

```
11.112.3.2 FactoryReset() void FactoryReset ( )
\textbf{11.112.3.3} \quad \textbf{GetBoardTemp()} \quad \texttt{unsigned int GetBoardTemp ()} \\
Gets the temperate of the mainboard in units of 0.1 °C.
11.112.3.4 GetCalibration() int GetCalibration (
               unsigned int channel )
11.112.3.5 GetCalibrationDecp() int GetCalibrationDecp (
               unsigned int channel )
\textbf{11.112.3.6} \quad \textbf{GetCalibrationMax()} \quad \texttt{int GetCalibrationMax} \quad \textbf{(}
               unsigned int channel )
11.112.3.7 GetCalibrationMin() int GetCalibrationMin (
               unsigned int channel )
11.112.3.8 GetCurrent() int GetCurrent (
               unsigned int channel )
11.112.3.9 GetD() int GetD (
               unsigned int device )
Gets the D-coefficient of the specified device.
11.112.3.10 GetDDecp() int GetDDecp (
               unsigned int device )
```

```
11.112.3.11 GetDevice() int GetDevice (
             unsigned int channel )
11.112.3.12 GetDeviceType() TcxDeviceTypeEnumNet GetDeviceType ( )
11.112.3.13 GetDevname() String ^ GetDevname (
             unsigned int device )
11.112.3.14 GetDMax() int GetDMax (
             unsigned int device )
11.112.3.15 GetDMin() int GetDMin (
             unsigned int device )
11.112.3.16 GetDuty() int GetDuty (
             unsigned int channel )
Gets the duty cycle of the heating element.
Parameters
 channel The channel number.
Returns
     The duty cycle in percent, the value of 320 \ast 64 corresponds to 100 %.
11.112.3.17 GetEnableHeaterLimit() bool GetEnableHeaterLimit (
             unsigned int device )
```

11.112.3.18 GetEnableThermocouple() bool GetEnableThermocouple (

unsigned int device)

```
\textbf{11.112.3.19} \quad \textbf{GetHasThermocouple()} \quad \texttt{bool GetHasThermocouple ()} \quad \texttt{bool GetHasThermocouple ()}
```

Gets weather the device supports a thermocouple.

```
11.112.3.20 GetHeaterLimit() int GetHeaterLimit ( unsigned int device )
```

Gets the temperate limit of the specified heater in units of $0.1\,^{\circ}$ C.

Gets the temperate of the specified heater in units of 0.1 °C.

```
11.112.3.22 Getl() int GetI (
unsigned int device)
```

Gets the I-coefficient of the specified device.

```
11.112.3.23 GetIDecp() int GetIDecp (
unsigned int device)
```

```
11.112.3.24 GetlMax() int GetlMax (
unsigned int device)
```

```
11.112.3.25 GetIMin() int GetIMin (
unsigned int device)
```

```
11.112.3.26 GetlOut() int GetIOut ( unsigned int channel)
```

Gets the current through the heating element.

channel The channel number.

Returns

The current in units of mA.

```
11.112.3.27 GetMaxHeaterPowerMultiwell() double GetMaxHeaterPowerMultiwell ( )
```

queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

```
11.112.3.28 GetMaxP() int GetMaxP (
unsigned int device)
```

Gets the maximum heater power of the specified device.

```
11.112.3.29 GetMaxpDecp() int GetMaxpDecp ( unsigned int device )
```

```
11.112.3.30 GetMaxpMax() int GetMaxpMax ( unsigned int device )
```

```
11.112.3.31 GetMaxpMin() int GetMaxpMin ( unsigned int device )
```

11.112.3.32 GetNumControlChannels() unsigned int GetNumControlChannels ()

Gets the number of channels the device can control/regulate.

11.112.3.33 GetNumDevices() unsigned int GetNumDevices ()

```
\textbf{11.112.3.34} \quad \textbf{GetNumMeasureChannels()} \quad \texttt{unsigned int GetNumMeasureChannels ()}
```

Gets the number of channels the device can measure.

```
11.112.3.35 GetOnOff() int GetOnOff (
          unsigned int channel )
```

Gets if the specified channel is on or off.

```
11.112.3.36 GetP() int GetP (
unsigned int device)
```

Gets the P-coefficient of the specified device.

```
11.112.3.37 GetPDecp() int GetPDecp ( unsigned int device )
```

```
11.112.3.38 GetPMax() int GetPMax (
unsigned int device)
```

```
11.112.3.39 GetPMin() int GetPMin (
unsigned int device)
```

```
11.112.3.40 GetPOut() int GetPOut (
unsigned int channel)
```

Gets the output power of the heating element.

Parameters

```
channel The channel number.
```

Returns

The resistance in units of mW.

```
11.112.3.41 GetPwrOut() int GetPwrOut (
             unsigned int channel )
11.112.3.42 GetPwrSet() int GetPwrSet (
             unsigned int channel )
11.112.3.43 GetRes1() int GetRes1 (
             unsigned int channel )
11.112.3.44 GetRes2() int GetRes2 (
             unsigned int channel )
11.112.3.45 GetResS() int GetResS (
             unsigned int channel )
11.112.3.46 GetResX() int GetResX (
             unsigned int channel )
11.112.3.47 GetROut() int GetROut (
             unsigned int channel )
Gets the resistance of the heating element.
Parameters
 channel The channel number.
```

Returns

The resistance in units of 0.1 Ohm.

```
11.112.3.48 GetSensorType() TcxSensorTypeEnumNet GetSensorType ( unsigned int device )
```

```
11.112.3.49 GetSetpoint() int GetSetpoint (
          unsigned int channel )
```

Gets the target temperate of specified channel in units of $0.1\,^{\circ}$ C.

```
11.112.3.53 GetThermocoupleCalibration() int GetThermocoupleCalibration ( unsigned int channel)
```

```
11.112.3.54 GetThermocoupleNanovoltPerKelvin() unsigned int GetThermocoupleNanovoltPerKelvin ( unsigned int channel)
```

Gets the proportional constant for the thermocouple.

Parameters

channel	Thermocouple channel number.
---------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

```
11.112.3.55 GetThermocoupleReferenceTemp() int GetThermocoupleReferenceTemp ( unsigned int channel)
```

Gets the voltage on the heating element.

Parameters

channel	The channel number.
---------	---------------------

Returns

The voltage in units of mV.

```
11.112.3.60 GetValue() int GetValue ( unsigned int channel )
```

Gets the temperate of the specified channel in units of 0.1 °C.

Gets the temperate of the specified channel in units of 0.01 ℃.

```
11.112.3.62 GetVolti() unsigned int GetVolti ( unsigned int channel )
```

```
11.112.3.63 SetCalibration() void SetCalibration (
             unsigned int channel,
             int calib )
11.112.3.64 SetD() void SetD (
             unsigned int device,
             int d\_coeff )
Sets the D-coefficient of the specified device.
11.112.3.65 SetDevice() void SetDevice (
             unsigned int channel,
             int device )
11.112.3.66 SetDeviceType() void SetDeviceType (
             TcxDeviceTypeEnumNet devicetype )
11.112.3.67 SetDevname() void SetDevname (
             unsigned int device,
             String^ Devicename )
11.112.3.68 SetEnableHeaterLimit() void SetEnableHeaterLimit (
             unsigned int device,
             bool enable )
11.112.3.69 SetEnableThermocouple() void SetEnableThermocouple (
             unsigned int device,
             bool enable )
11.112.3.70 SetHeaterLimit() void SetHeaterLimit (
             unsigned int device,
             int heater_limit )
```

```
11.112.3.71 Setl() void SetI (
         unsigned int device,
         int i_coeff )
```

Sets the I-coefficient of the specified device.

```
11.112.3.72 SetMaxHeaterPowerMultiwell() void SetMaxHeaterPowerMultiwell (

double MaxPowerWatt )
```

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

```
11.112.3.73 SetMaxP() void SetMaxP (
unsigned int device,
int maxp )
```

Sets the maximum heater power of the specified device.

```
11.112.3.74 SetOnOff() void SetOnOff (
     unsigned int channel,
     bool on )
```

Switches the specified channel on or off.

Parameters

```
channel The channel number.
```

```
11.112.3.75 SetP() void SetP ( unsigned int device, int p\_coeff )
```

Sets the P-coefficient of the specified device.

```
11.112.3.77 SetSetpoint() void SetSetpoint (
     unsigned int channel,
     int sp )
```

Sets the target temperate of specified channel in units of 0.1 ℃.

```
11.112.3.78 SetThermocoupleNanovoltPerKelvin() void SetThermocoupleNanovoltPerKelvin (
unsigned int channel,
unsigned int value)
```

Sets the proportional constant for the thermocouple.


Parameters

channel	Thermocouple channel number.
value	Proportinal constant in Nanovolt per Kelvin.

11.113 CTEERFunctionNet Class Reference

CTEERFunctionNet is the class to control the TEER device

Inheritance diagram for CTEERFunctionNet:

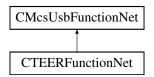

Public Member Functions

CTEERFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pTEERFunctionPointer
 — Container)

Initializes a new instance of the CTEERFunctionNet class.

- CTEERFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CTEERFunctionNet ()
- !CTEERFunctionNet ()
- uint32_t GetPeriod_us ()

gets the period of TEER stimulation in us

void SetPeriod_us (uint32_t period_us)

sets the period of TEER stimulation in us

• uint32 t GetAmplitude nA ()

gets TEER stimulation amplitude in nA

void SetAmplitude_nA (uint32_t Amplitude_nA)

sets TEER stimulation amplitude in nA

• TeerWaveformEnumNet GetWaveform ()

gets TEER stimulation waveform (sine/rect)

```
    void SetWaveform (TeerWaveformEnumNet Waveform)

      sets TEER stimulation waveform (sine/rect)

    TeerClampModeEnumNet GetClampMode ()

      gets TEER clamp mode (voltage/current)

    void SetClampMode (TeerClampModeEnumNet ClampMode)

     sets TEER clamp mode (voltage/current)

    void StartSampling (uint32_t NumberOfCycles)

      starts TEER stimulation (duration: n cycles) and samples during last cycle

    void StopSampling ()

     stops TEER stimulation and sampling

    uint32 t IsSamplingFinished ()

      returns false iff stimulation/sampling is going on, otherwise true

    void SetControllerParams (uint32 t P, uint32 t I, uint32 t D)

      sets PID controller parameters for voltage clamp mode

    void GetControllerParams ([System::Runtime::InteropServices::Out]uint32 t% P, [System::Runtime::←

  InteropServices::Out]uint32 t% I, [System::Runtime::InteropServices::Out]uint32 t% D)
      gets PID controller parameters for voltage clamp mode

    array< int32_t > ^ GetSampleBufferChunk (int Buffer_Length)

     private function to query max. 100 bytes of sample buffer; called internally

    array< int32_t > ^ GetSampleVoltageBuffer_uV (int Buffer_Length)

      returns voltage sample buffer (max. 500 values); unit: uV

    uint32 t GetMaxChunkSize Byte ()

     private function to be called internally only

    uint32 t GetBytesPerSample ()

     private function to be called internally only

    uint32 t GetNumberOfAvailableSamples ()

     private function to be called internally only

    void SetBufferIndex (uint32_t NewBufferIndex)

     pre-selects sample buffer to be tranferred by GetSampleVoltageBuffer uV()

    uint32 t GetAdapterCode ()

    uint32 t GetRotaryPositionCode ()

    void SetExternalLED (uint32 t NewState)

    void SetCurrentEnable (bool NewCurrentEnable)

      when disabled, no current will flow through chamber

    bool GetCurrentEnable ()

      when disabled, no current will flow through chamber

    int32_t GetUptimeSeconds ()

      returns time in seconds since device was powered up

    void StartInternalCalibration ()

      starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call
• bool IsInternalCalibrationFinished ()
      queries whether internal calibration has finished

    int GetDacZero ()

      returns DAC-offset (result of internal calibration); use to check for plausibility only

    void CancelInternalCalibration ()

     in case the internal calibration "hangs", this will cancel it
• void SetLiquidResistance (int32_t NewLiquidResistance_Ohm)

    int32_t GetLiquidResistance ()

    int GetScaleFactorU1 ()

      returns U1 scale factor times 10<sup>6</sup> (result of internal calibration)

    int GetScaleFactorU2 ()
```

```
returns U2 scale factor times 10^{\circ}6 (result of internal calibration)
```

int GetAdcOffsetU1 ()

returns ADC offset of U1 channel (result of internal calibration)

int GetAdcOffsetU2 ()

returns ADC offset of U2 channel (result of internal calibration)

int GetFrameErrorCounter ()

returns number of times (since bootup) sample memory got overwritten

• int GetSampleRate ()

returns sample rate in Hz

Additional Inherited Members

11.113.1 Detailed Description

CTEERFunctionNet is the class to control the TEER device

11.113.2 Constructor & Destructor Documentation

```
11.113.2.1 CTEERFunctionNet() [1/2] CTEERFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pTEERFunctionPointerContainer)
```

Initializes a new instance of the CTEERFunctionNet class.

```
11.113.2.2 CTEERFunctionNet() [2/2] CTEERFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.113.2.3 ~CTEERFunctionNet() virtual ~CTEERFunctionNet () [virtual]
```

```
11.113.2.4 "!CTEERFunctionNet() !CTEERFunctionNet ()
```

11.113.3 Member Function Documentation

```
11.113.3.1 CancelInternalCalibration() void CancelInternalCalibration ( )
in case the internal calibration "hangs", this will cancel it
11.113.3.2 GetAdapterCode() uint32_t GetAdapterCode ( )
Returns
11.113.3.3 GetAdcOffsetU1() int GetAdcOffsetU1 ( )
returns ADC offset of U1 channel (result of internal calibration)
Returns
11.113.3.4 GetAdcOffsetU2() int GetAdcOffsetU2 ( )
returns ADC offset of U2 channel (result of internal calibration)
Returns
11.113.3.5 GetAmplitude_nA() uint32_t GetAmplitude_nA ( )
gets TEER stimulation amplitude in nA
Returns
     current stimulation amplitude in nA
11.113.3.6 GetBytesPerSample() uint32_t GetBytesPerSample ()
private function to be called internally only
Returns
```

```
\textbf{11.113.3.7} \quad \textbf{GetClampMode()} \quad \texttt{TeerClampModeEnumNet GetClampMode ( )}
```

gets TEER clamp mode (voltage/current)

Returns

current TEER clamp mode

```
11.113.3.8 GetControllerParams() void GetControllerParams (
```

```
[System::Runtime::InteropServices::Out] uint32_t% P,
[System::Runtime::InteropServices::Out] uint32_t% I,
[System::Runtime::InteropServices::Out] uint32_t% D)
```

gets PID controller parameters for voltage clamp mode

Parameters

Р	
1	
D	

11.113.3.9 GetCurrentEnable() bool GetCurrentEnable ()

when disabled, no current will flow through chamber

Returns

```
11.113.3.10 GetDacZero() int GetDacZero ( )
```

returns DAC-offset (result of internal calibration); use to check for plausibility only

Returns

11.113.3.11 GetFrameErrorCounter() int GetFrameErrorCounter ()

returns number of times (since bootup) sample memory got overwritten

Returns

```
11.113.3.12 GetLiquidResistance() int32_t GetLiquidResistance ( )
Returns
11.113.3.13 GetMaxChunkSize_Byte() uint32_t GetMaxChunkSize_Byte ( )
private function to be called internally only
Returns
11.113.3.14 GetNumberOfAvailableSamples() uint32_t GetNumberOfAvailableSamples ( )
private function to be called internally only
Returns
11.113.3.15 GetPeriod_us() uint32_t GetPeriod_us ( )
gets the period of TEER stimulation in us
Returns
11.113.3.16 GetRotaryPositionCode() uint32_t GetRotaryPositionCode ( )
Returns
11.113.3.17 GetSampleBufferChunk() array<int32_t> ^{\land} GetSampleBufferChunk (
              int Buffer_Length )
private function to query max. 100 bytes of sample buffer; called internally
```

Pa	ra	m	ρį	ŀΔ	re

Buffer_Length	The maximal length of Buffer.
---------------	-------------------------------

Returns

11.113.3.18 GetSampleRate() int GetSampleRate ()

returns sample rate in Hz

Returns

11.113.3.19 GetSampleVoltageBuffer_uV() array<int32_t> $^{\land}$ GetSampleVoltageBuffer_uV (int Buffer_Length)

returns voltage sample buffer (max. 500 values); unit: uV

Parameters

Buffer Lenath	The maximal length of Buffer.

Returns

11.113.3.20 GetScaleFactorU1() int GetScaleFactorU1 ()

returns U1 scale factor times $10^{\circ}6$ (result of internal calibration)

Returns

```
11.113.3.21 GetScaleFactorU2() int GetScaleFactorU2 ( )
returns U2 scale factor times 10<sup>6</sup> (result of internal calibration)
Returns
11.113.3.22 GetUptimeSeconds() int32_t GetUptimeSeconds ( )
returns time in seconds since device was powered up
Returns
11.113.3.23 GetWaveform() TeerWaveformEnumNet GetWaveform ( )
gets TEER stimulation waveform (sine/rect)
Returns
     waveform enum
11.113.3.24 IsInternalCalibrationFinished() bool IsInternalCalibrationFinished ( )
queries whether internal calibration has finished
Returns
11.113.3.25 IsSamplingFinished() uint32_t IsSamplingFinished ( )
returns false iff stimulation/sampling is going on, otherwise true
Returns
11.113.3.26 SetAmplitude_nA() void SetAmplitude_nA (
              uint32_t Amplitude_nA )
```

sets TEER stimulation amplitude in nA

Amplitude_nA	new stimulation amplitude in nA	
--------------	---------------------------------	--

```
11.113.3.27 SetBufferIndex() void SetBufferIndex ( uint32_t NewBufferIndex )
```

pre-selects sample buffer to be tranferred by GetSampleVoltageBuffer_uV()

Parameters

NewBufferIndex 0 - chamber voltage; 1 - compliance voltage

```
11.113.3.28 SetClampMode() void SetClampMode (
TeerClampModeEnumNet ClampMode)
```

sets TEER clamp mode (voltage/current)

Parameters

ClampMode	new TEER clamp mode
-----------	---------------------

sets PID controller parameters for voltage clamp mode

Parameters

Р	
1	
D	

```
11.113.3.30 SetCurrentEnable() void SetCurrentEnable ( bool NewCurrentEnable )
```

when disabled, no current will flow through chamber

NewCurrentEnable

```
11.113.3.31 SetExternalLED() void SetExternalLED ( uint32_t NewState )
```

Parameters

NewState

```
11.113.3.32 SetLiquidResistance() void SetLiquidResistance ( int32_t NewLiquidResistance_Ohm )
```

Parameters

NewLiquidResistance_Ohm

```
11.113.3.33 SetPeriod_us() void SetPeriod_us ( uint32_t period_us )
```

sets the period of TEER stimulation in us

Parameters

period_us

```
11.113.3.34 SetWaveform() void SetWaveform (

TeerWaveformEnumNet Waveform)
```

sets TEER stimulation waveform (sine/rect)

Parameters

Waveform waveform enum

$\textbf{11.113.3.35} \quad \textbf{StartInternalCalibration()} \quad \texttt{void StartInternalCalibration ()}$

starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call

```
11.113.3.36 StartSampling() void StartSampling ( uint32_t NumberOfCycles )
```

starts TEER stimulation (duration: n cycles) and samples during last cycle

Parameters


NumberOfCycles	number of cycles (sine or rect) to output (0 - loop forever)
----------------	--

11.113.3.37 StopSampling() void StopSampling ()

stops TEER stimulation and sampling

11.114 CTEERMachineDeviceNet Class Reference

Inheritance diagram for CTEERMachineDeviceNet:

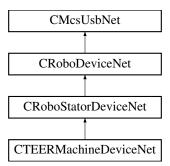

Public Member Functions

- CTEERMachineDeviceNet ()
- ∼CTEERMachineDeviceNet ()

Properties

• CTEERFunctionNet[^] TEERFunctionNet [get]

Additional Inherited Members

11.114.1 Constructor & Destructor Documentation

11.114.1.1 CTEERMachineDeviceNet() CTEERMachineDeviceNet ()

11.114.1.2 \sim CTEERMachineDeviceNet() \sim CTEERMachineDeviceNet ()

11.114.2 Property Documentation

11.114.2.1 TEERFunctionNet CTEERFunctionNet^ TEERFunctionNet [get]

11.115 CUsbExceptionNet Class Reference

Exception class that is thrown in case of an USB error.

Inheritance diagram for CUsbExceptionNet:



Public Member Functions

- CUsbExceptionNet (uint32_t status)

 Constructor of a CUsbException.
- CUsbExceptionNet (uint32_t status, String^ message)

Properties

• uint32_t Status [get]

11.115.1 Detailed Description

Exception class that is thrown in case of an USB error.

11.115.2 Constructor & Destructor Documentation

```
11.115.2.1 CUsbExceptionNet() [1/2] CUsbExceptionNet ( uint32_t status )
```

Constructor of a CUsbException.

status the status number

11.115.3 Property Documentation

```
11.115.3.1 Status uint32_t Status [get]
```

11.116 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference

Public Member Functions

CVoltageRangeInfoNet (int vr, String[^] vrString)

Public Attributes

- int VoltageRangeInMicroVolt
- String ^ VoltageRangeDisplayStringMilliVolt

11.116.1 Constructor & Destructor Documentation


11.116.2 Member Data Documentation

$\textbf{11.116.2.1} \quad \textbf{VoltageRangeDisplayStringMilliVolt} \quad \texttt{String} \; \wedge \; \texttt{VoltageRangeDisplayStringMilliVolt}$

11.116.2.2 VoltageRangeInMicroVolt int VoltageRangeInMicroVolt

11.117 CW2100_FunctionNet Class Reference

Inheritance diagram for CW2100 FunctionNet:

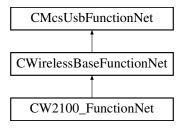

Classes

struct AudioChannelsNet

Public Member Functions

- CW2100_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] w2100_Function
 — PointerContainer)
- CW2100_FunctionNet (CMcsUsbNet[∧] mcsusb)
- array< HeadStageIDType[^]> [^] GetAvailableHeadstages (unsigned int max_length)
- void SelectHeadstage (unsigned int IDorEntry, int TimeSlotNr)
- void DeselectHeadstage (int TimeSlotNr)
- void DeselectAllHeadstages ()
- HeadStageIDTypeState ^ GetSelectedHeadstageState (int TimeSlotNr)
- BatteryState ^ GetBatteryState (int TimeSlotNr)
- System::String \(^\) GetUserDefinedName (unsigned short ID)
- System::String ^ GetUserDefinedNameFromSelectedHS (int TimeSlotNr)
- System::String \(^\) GetUserDefinedNameCache (unsigned short ID)
- uint32_t GetUserDefinedNameCache (unsigned short ID, [System::Runtime::InteropServices::Out]System

 ::String^% Name)
- W2100_StimulusParametersNet ^ GetStiumlusParameters (unsigned short ID)
- W2100 StimulusParametersNet ^ GetStimulusParametersFromSelectedHS (int TimeSlotNr)
- W2100 StimulusParametersNet \(^\) GetStimulusParametersCache (unsigned int typeValue)
- uint32_t GetStimulusParametersCache (unsigned int typeValue, [System::Runtime::InteropServices::
 — Out]W2100_StimulusParametersNet^% StimulusParameters)
- void SetSelectedChannels (array< BYTE >^ channels, int TimeSlotNr)
- array< BYTE > ^ GetSelectedChannels (int TimeSlotNr)
- void SetMultiHeadstageMode (bool Mode)
- bool GetMultiHeadstageMode ()
- void SetHeadstageSamplingActive (bool Active, int TimeSlotNr)
- bool GetHeadstageSamplingActive (int TimeSlotNr)
- void SetHeadstageToSleep (unsigned int Sleep16ms, int TimeSlotNr)
- void SetHeadstageOnOff (unsigned short On, int TimeSlotNr)
- unsigned short GetHeadstageOnOff (int TimeSlotNr)
- unsigned int GetAnalogOutChannel ([System::Runtime::InteropServices::Out]int % automatic, unsigned short index)
- void SetAnalogOutChannel (int automatic, unsigned short index, unsigned int Channel)

- array< unsigned int > ^ GetAnalogOutFilter ([System::Runtime::InteropServices::Out]int % automatic)
- void SetAnalogOutFilter (int automatic, array< unsigned int >^ Coeffs)
- AnalogOut_DAC_Range_EnumNet GetDacRange ()
- void SetDacRange (AnalogOut_DAC_Range_EnumNet range)
- CFilterPropertyNet ^ GetFilterProperty (W2100DacqGroupChannelEnumNet GroupID, unsigned int index)
- array < CFilterPropertyNet^> ^ GetFilterProperties (W2100DacqGroupChannelEnumNet GroupID)
- void SetAccelGyroEnabled (W2100_Accel_Gyro_Select_EnumNet enable, int TimeSlotNr)
- W2100 Accel Gyro Select EnumNet GetAccelGyroEnabled (int TimeSlotNr)
- void SetAccelGyroDesiredRate (int rate, int TimeSlotNr)
- int GetAccelGyroDesiredRate (int TimeSlotNr)
- int GetAccelGyroCurrentRate (int TimeSlotNr)
- void SetAccelRange (int range, int TimeSlotNr)
- int GetAccelRange (int TimeSlotNr)
- void SetGyroRange (int range, int TimeSlotNr)
- int GetGyroRange (int TimeSlotNr)
- void SetAudioChannels (array< AudioChannelsNet[^]>[^] channels)
- array< AudioChannelsNet[^]> [^] GetAudioChannels ()
- unsigned int GetPicFirmwareType (int TimeSlotNr)
- unsigned int GetFPGAFirmwareType (int TimeSlotNr)

Static Public Member Functions

- static void ClearUserDefinedNameCache ()
- static void ClearUserDefinedNameCache (unsigned short ID)
- static void ClearStimulusParametersCache ()
- static void ClearStimulusParametersCache (unsigned short ID)

Properties

- CW2100_StimulatorFunctionNet^ Stimulator [get]
- CPulseGeneratorFunctionNet[^] PulseGenerator [get]

Additional Inherited Members

11.117.1 Constructor & Destructor Documentation

11.117.2 Member Function Documentation

```
11.117.2.1 ClearStimulusParametersCache() [1/2] static void ClearStimulusParametersCache ( )
[static]
11.117.2.2 ClearStimulusParametersCache() [2/2] static void ClearStimulusParametersCache (
              unsigned short ID ) [static]
11.117.2.3 ClearUserDefinedNameCache() [1/2] static void ClearUserDefinedNameCache ( ) [static]
\textbf{11.117.2.4} \quad \textbf{ClearUserDefinedNameCache() [2/2]} \quad \texttt{static void ClearUserDefinedNameCache ()} \\
              unsigned short ID ) [static]
11.117.2.5 DeselectAllHeadstages() void DeselectAllHeadstages ()
11.117.2.6 DeselectHeadstage() void DeselectHeadstage (
              int TimeSlotNr )
\textbf{11.117.2.7} \quad \textbf{GetAccelGyroCurrentRate()} \quad \texttt{int GetAccelGyroCurrentRate ()}
              int TimeSlotNr )
11.117.2.8 GetAccelGyroDesiredRate() int GetAccelGyroDesiredRate (
              int TimeSlotNr )
11.117.2.9 GetAccelGyroEnabled() W2100_Accel_Gyro_Select_EnumNet GetAccelGyroEnabled (
              int TimeSlotNr )
```

```
11.117.2.10 GetAccelRange() int GetAccelRange (
             int TimeSlotNr )
11.117.2.11 GetAnalogOutChannel() unsigned int GetAnalogOutChannel (
             [System::Runtime::InteropServices::Out] int % automatic,
             unsigned short index )
11.117.2.12 GetAnalogOutFilter() array<unsigned int> ^ GetAnalogOutFilter (
             [System::Runtime::InteropServices::Out] int % automatic )
11.117.2.13 GetAudioChannels() array<AudioChannelsNet^> ^ GetAudioChannels ( )
11.117.2.14 GetAvailableHeadstages() array<HeadStageIDType^> ^ GetAvailableHeadstages (
             unsigned int max_length )
11.117.2.15 GetBatteryState() BatteryState ^ GetBatteryState (
             int TimeSlotNr )
11.117.2.16 GetDacRange() AnalogOut_DAC_Range_EnumNet GetDacRange ( )
11.117.2.17 GetFilterProperties() array<CFilterPropertyNet^> ^ GetFilterProperties (
             W2100DacqGroupChannelEnumNet GroupID )
11.117.2.18 GetFilterProperty() CFilterPropertyNet ^ GetFilterProperty (
             W2100DacqGroupChannelEnumNet GroupID,
             unsigned int index )
11.117.2.19 GetFPGAFirmwareType() unsigned int GetFPGAFirmwareType (
             int TimeSlotNr )
```

```
\textbf{11.117.2.20} \quad \textbf{GetGyroRange()} \quad \texttt{int GetGyroRange ()}
                                                                    int TimeSlotNr )
11.117.2.21 GetHeadstageOnOff() unsigned short GetHeadstageOnOff (
                                                                   int TimeSlotNr )
11.117.2.22 GetHeadstageSamplingActive() bool GetHeadstageSamplingActive (
                                                                    int TimeSlotNr )
11.117.2.23 GetMultiHeadstageMode() bool GetMultiHeadstageMode ( )
11.117.2.24 GetPicFirmwareType() unsigned int GetPicFirmwareType (
                                                                   int TimeSlotNr )
11.117.2.25 GetSelectedChannels() array<br/>
BYTE> ^{\wedge} GetSelectedChannels (
                                                                    int TimeSlotNr )
\textbf{11.117.2.26} \quad \textbf{GetSelectedHeadstageState()} \quad \texttt{HeadStageIDTypeState} \quad ^{\land} \quad \texttt{GetSelectedHeadstageState} \quad \texttt{(}
                                                                    int TimeSlotNr )
11.117.2.27 GetStimulusParametersCache() [1/2] W2100_StimulusParametersNet ^ GetStimulus↔
ParametersCache (
                                                                   unsigned int typeValue )
\textbf{11.117.2.28} \quad \textbf{GetStimulusParametersCache() [2/2]} \quad \texttt{uint32\_t GetStimulusParametersCache ()} \quad \textbf{(2/2)} \quad \textbf{(2/2)
                                                                  unsigned int typeValue,
                                                                    [System::Runtime::InteropServices::Out] W2100_StimulusParametersNet^% Stimulus↔
Parameters )
```


```
11.117.2.29 GetStimulusParametersFromSelectedHS() w2100_StimulusParametersNet ^ GetStimulus↔
ParametersFromSelectedHS (
             int TimeSlotNr )
11.117.2.30 GetStiumlusParameters() W2100_StimulusParametersNet ^ GetStiumlusParameters (
             unsigned short ID )
\textbf{11.117.2.31} \quad \textbf{GetUserDefinedName()} \quad \texttt{System::String} \; \wedge \; \texttt{GetUserDefinedName} \; \; (
             unsigned short ID )
11.117.2.32 GetUserDefinedNameCache() [1/2] System::String ^ GetUserDefinedNameCache (
             unsigned short ID )
11.117.2.33 GetUserDefinedNameCache() [2/2] uint32_t GetUserDefinedNameCache (
             unsigned short ID,
             [System::Runtime::InteropServices::Out] System::String^{\land}% Name )
11.117.2.34 GetUserDefinedNameFromSelectedHS() System::String ^ GetUserDefinedNameFrom←
SelectedHS (
             int TimeSlotNr )
11.117.2.35 SelectHeadstage() void SelectHeadstage (
             unsigned int IDorEntry,
             int TimeSlotNr )
11.117.2.36 SetAccelGyroDesiredRate() void SetAccelGyroDesiredRate (
             int rate,
             int TimeSlotNr )
11.117.2.37 SetAccelGyroEnabled() void SetAccelGyroEnabled (
             W2100_Accel_Gyro_Select_EnumNet enable,
             int TimeSlotNr )
```

```
11.117.2.38 SetAccelRange() void SetAccelRange (
             int range,
             int TimeSlotNr )
11.117.2.39 SetAnalogOutChannel() void SetAnalogOutChannel (
            int automatic,
             unsigned short index,
             unsigned int Channel )
11.117.2.40 SetAnalogOutFilter() void SetAnalogOutFilter (
             int automatic,
             array< unsigned int >^{\land} Coeffs )
11.117.2.41 SetAudioChannels() void SetAudioChannels (
             array< AudioChannelsNet^>^ channels )
11.117.2.42 SetDacRange() void SetDacRange (
             AnalogOut_DAC_Range_EnumNet range )
11.117.2.43 SetGyroRange() void SetGyroRange (
             int range,
             int TimeSlotNr )
11.117.2.44 SetHeadstageOnOff() void SetHeadstageOnOff (
             unsigned short On,
             int TimeSlotNr )
11.117.2.45 SetHeadstageSamplingActive() void SetHeadstageSamplingActive (
            bool Active,
             int TimeSlotNr )
```

11.118 CW2100_StimulatorFunctionNet Class Reference

11.117.3.2 Stimulator CW2100_StimulatorFunctionNet^ Stimulator [get]

Inheritance diagram for CW2100_StimulatorFunctionNet:

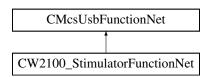

Public Member Functions

- CW2100 StimulatorFunctionNet (CMcsUsbNet[^] mcsusb)
- void SendStart (uint32_t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

void SendStop (uint32_t triggermap)

Stop some or all triggers of the STG.

- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareData (int channel, array< int32_t > ^ amplitude, array< uint64_t > ^ duration, STG_DestinationEnumNet destType, uint32_t repeat)
- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareDataSync (int channel, array< int32_t >^ amplitude, array< uint32_t >^ Sync, array< uint64_t >^ duration, STG_DestinationEnumNet destType, uint32_t repeat)
- void SendPreparedData (int channel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled, STG_DestinationEnumNet destType)
- void ClearChannelData (int channel)

Delete a Stimulus Pattern from STG memory

• int GetDACResolution ()

Gets number of bits of the DAC resolution.

int GetTimeResolutionInNanoSeconds ()

Gets number of bits of the DAC resolution.

int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

• int GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

int GetCurrentRangeInNanoAmp (uint32 t channel)

Gets the Current Range of the specified channel in Nanoamps.

int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

- uint32 t GetNumberOfAnalogChannels ()
- uint32_t GetNumberOfSyncoutChannels ()
- uint32 t GetNumberOfTriggerInputs ()
- void SelectTimeSlot (int TimeSlotNr)
- int GetTimeSlot ()
- uint32_t GetStimulationPatternMemory ()
- uint32 t GetBoostPreTime ()
- void SetBoostAlwaysOnMode (uint32_t alwayson_mode)
- uint32_t GetBoostAlwaysOnMode ()
- void SetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, W2100DigitalSourceEnumNet digstream_source, int bitnumber_offset)
- void GetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% digstream_source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)
- void SetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, DigitalStimulatorTriggerSlopeEnumNet slope)
- DigitalStimulatorTriggerSlopeEnumNet GetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulator
 — TriggerEventEnumNet trigger_event, int trigger_number)
- void StartPoll ()
- void StopPoll ()

Static Public Attributes

- static const uint32 t BOOST BIT = (1 << 0)
- static const uint32_t GND_SWITCH_BIT = (1 << 1)
- static const uint32_t SYNC_BIT0 = (1 << 2)
- static const uint32_t SYNC_BIT1 = (1 << 3)

Events

OnStgPollStatus[^] PollStatusEvent

Additional Inherited Members

11.118.1 Constructor & Destructor Documentation

```
11.118.1.1 CW2100_StimulatorFunctionNet() CW2100_StimulatorFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.118.2 Member Function Documentation

Delete a Stimulus Pattern from STG memory

Parameters

channel specifies the channel to clear.

```
11.118.2.2 GetBoostAlwaysOnMode() uint32_t GetBoostAlwaysOnMode ( )
```

```
11.118.2.3 GetBoostPreTime() uint32_t GetBoostPreTime ( )
```

```
11.118.2.4 GetCurrentRangeInNanoAmp() int GetCurrentRangeInNanoAmp ( uint32_t channel)
```

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel Channel which is queried.

Returns

The Current Range of the specified channel in Nanoamps.

11.118.2.5 GetCurrentResolutionInNanoAmp() int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.118.2.6 GetDACResolution() int GetDACResolution ()

Gets number of bits of the DAC resolution.

Returns

source.

The DAC resolution in bits.

```
11.118.2.7 GetDigitalStimulatorTrigger() void GetDigitalStimulatorTrigger (
```

```
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
[System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% digstream_←
[System::Runtime::InteropServices::Out] int% bitnumber_offset )
```

11.118.2.8 GetDigitalStimulatorTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetDigital←

```
\textbf{11.118.2.9} \quad \textbf{GetNumberOfAnalogChannels()} \quad \texttt{uint32\_t GetNumberOfAnalogChannels ()} \\
11.118.2.10 GetNumberOfSyncoutChannels() uint32_t GetNumberOfSyncoutChannels ( )
11.118.2.11 GetNumberOfTriggerInputs() uint32_t GetNumberOfTriggerInputs ( )
11.118.2.12 GetStimulationPatternMemory() uint32_t GetStimulationPatternMemory ( )
11.118.2.13 GetTimeResolutionInNanoSeconds() int GetTimeResolutionInNanoSeconds ( )
Gets number of bits of the DAC resolution.
Returns
     The time resolution in ns.
11.118.2.14 GetTimeSlot() int GetTimeSlot ( )
11.118.2.15 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (
              uint32_t channel )
Gets the Voltage Range of the specified channel in Microvolts.
Parameters
 channel
           Channel which is queried.
```

Generated by Doxygen

The Voltage Range of the specified channel in Microvolts.

Returns

11.118.2.16 GetVoltageResolutionInMicroVolt() int GetVoltageResolutionInMicroVolt (uint32_t *channel*)

Gets the Voltage Resolution of the specified channel in Microvolts.

channel Chai	nnel which is queried.
--------------	------------------------

Returns

The Voltage Resolution of the specified channel in Microvolts.

```
11.118.2.17 PrepareData() CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ Prepare↔
Data (
              int channel,
              array< int32_t >^{\land} amplitude,
              array< uint64_t >^{\land} duration,
              STG_DestinationEnumNet destType,
              uint32_t repeat )
11.118.2.18 PrepareDataSync() CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ Prepare←
DataSync (
             int channel,
             array< int32_t >^{\wedge} amplitude,
              array< uint32_t >^{\land} Sync,
              array< uint64_t >^{\wedge} duration,
              STG_DestinationEnumNet destType,
              uint32_t repeat )
11.118.2.19 SelectTimeSlot() void SelectTimeSlot (
              int TimeSlotNr )
11.118.2.20 SendPreparedData() void SendPreparedData (
             int channel,
              CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^ device_data_and_unrolled,
              {\tt STG\_DestinationEnumNet}\ \textit{destType}\ )
11.118.2.21 SendStart() void SendStart (
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

uint32_t triggermap)

triggermap	A bitmap of triggers which will be started.
------------	---

```
11.118.2.22 SendStop() void SendStop ( uint32_t triggermap)
```

Stop some or all triggers of the STG.

Parameters

triggermap	A bitmap of triggers which will be stopped.
------------	---

```
11.118.2.23 SetBoostAlwaysOnMode() void SetBoostAlwaysOnMode ( uint32_t alwayson_mode )
```

```
11.118.2.24 SetDigitalStimulatorTrigger() void SetDigitalStimulatorTrigger (
```

```
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
W2100DigitalSourceEnumNet digstream_source,
int bitnumber_offset )
```

$\textbf{11.118.2.25} \quad \textbf{SetDigitalStimulatorTriggerSlope()} \quad \texttt{void SetDigitalStimulatorTriggerSlope} \quad \textbf{(}$

```
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
DigitalStimulatorTriggerSlopeEnumNet slope )
```

```
11.118.2.26 StartPoll() void StartPoll ()
```

$\textbf{11.118.2.27} \quad \textbf{StopPoll()} \quad \texttt{void StopPoll ()}$

11.118.3 Member Data Documentation

```
11.118.3.1 BOOST_BIT const uint32_t BOOST_BIT = (1 << 0) [static]
```

```
\textbf{11.118.3.2} \quad \textbf{GND\_SWITCH\_BIT} \quad \texttt{const uint32\_t GND\_SWITCH\_BIT} = (1 << 1) \quad \texttt{[static]}
```

```
11.118.3.3 SYNC_BIT0 const uint32_t SYNC_BIT0 = (1 << 2) [static]
```


```
11.118.3.4 SYNC_BIT1 const uint32_t SYNC_BIT1 = (1 << 3) [static]
```

11.118.4 Event Documentation

11.118.4.1 PollStatusEvent OnStgPollStatus^ PollStatusEvent

11.119 CW2100DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CW2100DacqGroupChannelSelectionNet:

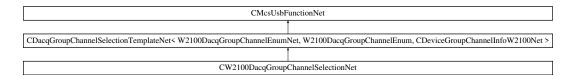

Public Member Functions

CW2100DacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members


11.119.1 Constructor & Destructor Documentation

$\textbf{11.119.1.1} \quad \textbf{CW2100DacqGroupChannelSelectionNet()} \quad \texttt{CW2100DacqGroupChannelSelectionNet} \quad (\\ & \texttt{CMcsUsbNet} \land \ \textit{mcsusb} \)$

11.120 CWarnerUssingDeviceNet Class Reference

CWarnerUssingDeviceNet is the class to control the Ussing device

Inheritance diagram for CWarnerUssingDeviceNet:

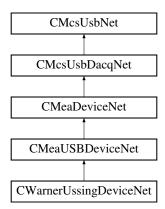

Public Member Functions

- CWarnerUssingDeviceNet ()
 - Initializes a new instance of the CWarnerUssingDeviceNet class.
- virtual ~CWarnerUssingDeviceNet ()
- !CWarnerUssingDeviceNet ()

Properties

CWarnerUssingFunctionNet[^] WarnerUssingFunction [get]

Additional Inherited Members

11.120.1 Detailed Description

CWarnerUssingDeviceNet is the class to control the Ussing device

11.120.2 Constructor & Destructor Documentation

11.120.2.1 CWarnerUssingDeviceNet() CWarnerUssingDeviceNet ()

Initializes a new instance of the CWarnerUssingDeviceNet class.

$\textbf{11.120.2.2} \quad \sim \textbf{CWarnerUssingDeviceNet()} \quad \text{virtual} \quad \sim \texttt{CWarnerUssingDeviceNet} \quad \textbf{()} \quad \text{[virtual]}$

11.120.2.3 "!CWarnerUssingDeviceNet() !CWarnerUssingDeviceNet ()


11.120.3 Property Documentation

11.120.3.1 WarnerUssingFunction CWarnerUssingFunctionNet^ WarnerUssingFunction [get]

11.121 CWarnerUssingFunctionNet Class Reference

CWarnerUssingFunctionNet is the class to control the Ussing device

Inheritance diagram for CWarnerUssingFunctionNet:

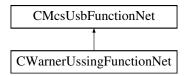

Public Member Functions

CWarnerUssingFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pWarner
 — UssingFunctionPointerContainer)

Initializes a new instance of the CWarnerUssingFunctionNet class.

- CWarnerUssingFunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CWarnerUssingFunctionNet ()
- !CWarnerUssingFunctionNet ()
- int32_t GetChannelsCountOfChamber (int32_t ChamberId)

gets number of channels in datastream from chamber amp with given index

int32_t GetNumberOfHardwareSlotsForChambers ()

gets number of physical hardware slots for chambers amps

int32 t GetNumberOfAvailableChambers ()

gets number of actually connected chamber amps

bool IsChamberAvailable (int32_t ChamberId)

checks whether chamber amp is connected to slot

 void SetPulse (int32_t ChamberId, UssingClampModeEnumNet StgMode, int32_t NumberOfRepetitions, array< int >^ Amplitudes, array< int >^ Durations)

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

void SetVoltageClampControllerParam_P (int32_t ChamberId, uint32_t P)

sets P value of PID controller;

• void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)

sets I value of PID controller;

void SetVoltageClampControllerParam_D (int32_t ChamberId, uint32_t D)

sets D value of PID controller;

uint32_t GetVoltageClampControllerParam_P (int32_t ChamberId)

gets P value of PID controller;

```
    uint32 t GetVoltageClampControllerParam I (int32 t ChamberId)

     gets I value of PID controller;

    uint32_t GetVoltageClampControllerParam_D (int32_t ChamberId)

     gets D value of PID controller:

    void SetClampMode (int32 t Chamberld, UssingClampModeEnumNet NewClampMode)

     sets clamp mode (voltage, current or open clamp)

    UssingClampModeEnumNet GetClampMode (int32 t ChamberId)

     gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)

    bool IsInternalCalibrationFinished (int32 t ChamberId)

      when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

    int32 t GetU1Offset (int32 t ChamberId)

    int32 t GetU2Offset (int32 t ChamberId)

    int32_t GetU1Reference (int32_t Chamberld)

    int32_t GetU2Reference (int32_t ChamberId)

    int32 t GetDacZero (int32 t ChamberId)

    void SetHighCurrentMode (int32 t ChamberId)

    void SetLowCurrentMode (int32 t ChamberId)

    bool IsHighCurrentMode (int32_t ChamberId)

    uint32_t GetLowCurrentRange (int32_t Chamberld)

    uint32 t GetHighCurrentRange (int32 t Chamberld)

     unit: nA

    uint32_t GetDacPampsPerDigitLowCurrentRange (int32_t ChamberId)

     unit: pA/digit

    uint32 t GetDacNampsPerDigitHighCurrentRange (int32 t ChamberId)

     unit: nA/digit

    uint32_t GetUnitsPerDigit (int32_t Chamberld, int32_t Channelld)

     gets amps/volts per digit for specified chamber and channel

    int32 t GetUnitExponent (int32 t Chamberld, int32 t Channelld)

    UssingUnitEnumNet GetUnitName (int32_t ChamberId, int32_t ChannelId)

     gets the channel's unit name

    String \(^\) GetUnitDescription (int32 t Chamberld, int32 t Channelld)

    array< int > ^ GetAvailableChambers ()

     returns array with (zero-based) Chamberlds of all available chambers

    int32_t GetUptimeSeconds (int32_t ChamberId)

    void SetIdleModeOffset (int32 t ChamberId, UssingClampModeEnumNet ClampMode, int32 t NewIdle

  Offset)
      sets the offset (voltage or current) that will be applied when clamping is DISABLED

    int32 t GetIdleModeOffset (int32 t ChamberId, UssingClampModeEnumNet ClampMode)

     gets the offset (voltage or current) that will be applied when clamping is DISABLED

    void SetEnablePulse (int32_t Chamberld, UssingClampModeEnumNet ClampMode, bool Enable)

     enable pulse of given chamber and mode (voltage/current clamp) of this chamber

    bool IsPulseEnabled (int32 t ChamberId, UssingClampModeEnumNet ClampMode)

     returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED

    void SetLiquidResistance (int32 t Chamberld, int32 t NewLiquidResistance Ohm)

    int32_t GetLiquidResistance (int32_t ChamberId)

    int32 t GetComplianceVoltage (int32 t ChamberId)

     returns compliance voltage threshold in digits; when Uc is above, current source is overloaded

    bool CompensateElectrodeOffset (int32_t ChamberId)
```

blocking call to compensate electrode offset of one chamber; returns true when successful

blocking call that waits for chamber boot-up calibration to complete

blocking call that waits for ALL chambers' boot-up calibration to complete

Generated by Doxygen

bool WaitForChamber (int32_t ChamberId)

bool WaitForAllChambers ()

Additional Inherited Members

11.121.1 Detailed Description

CWarnerUssingFunctionNet is the class to control the Ussing device

11.121.2 Constructor & Destructor Documentation

```
11.121.2.1 CWarnerUssingFunctionNet() [1/2] CWarnerUssingFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pWarnerUssingFunctionPointerContainer)
```

Initializes a new instance of the CWarnerUssingFunctionNet class.

```
11.121.2.2 CWarnerUssingFunctionNet() [2/2] CWarnerUssingFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.121.2.3 ~CWarnerUssingFunctionNet() virtual ~CWarnerUssingFunctionNet () [virtual]
```

```
11.121.2.4 "!CWarnerUssingFunctionNet() !CWarnerUssingFunctionNet ( )
```

11.121.3 Member Function Documentation

```
11.121.3.1 CompensateElectrodeOffset() bool CompensateElectrodeOffset ( int32_t ChamberId )
```

blocking call to compensate electrode offset of one chamber; returns true when successful

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.2 GetAvailableChambers() array<int> $^{\wedge}$ GetAvailableChambers ()

returns array with (zero-based) Chamberlds of all available chambers

gets number of channels in datastream from chamber amp with given index

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

return value of zero means that amp is not placed

11.121.3.4 GetClampMode() UssingClampModeEnumNet GetClampMode (int32_t ChamberId)

gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

returns compliance voltage threshold in digits; when Uc is above, current source is overloaded

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

unit: nA/digit

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.7 GetDacPampsPerDigitLowCurrentRange() uint32_t GetDacPampsPerDigitLowCurrentRange (int32_t ChamberId)

unit: pA/digit

Parameters

01 1	
Chamber⊸	index of hardware chamber slot (zero-based)
ld	

Returns

• diagnostic function only - ; gets real zero value of DAC in digits (0 -> neg. current; 32767 -> near zero; 65535 -> pos. current)

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

unit: nA

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

gets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

Chamberld	index of hardware chamber slot (zero-based)
ClampMode	voltage or current clamp stimulation

Returns

unit: nA or mV

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

unit: nA

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

 $\textbf{11.121.3.13} \quad \textbf{GetNumberOfAvailableChambers()} \quad \texttt{int32_t} \quad \texttt{GetNumberOfAvailableChambers} \quad \textbf{()}$

gets number of actually connected chamber amps

Returns

11.121.3.14 GetNumberOfHardwareSlotsForChambers() int32_t GetNumberOfHardwareSlotsFor← Chambers ()

gets number of physical hardware slots for chambers amps

Returns

· diagnostic function only -

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

· diagnostic function only -

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

· diagnostic function only -

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

· diagnostic function only -

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

Parameters

Chamber← Id	index of hardware chamber slot (zero-based)
Channelld	index of channel (zero-based)

Returns

Parameters

Chamber← Id	index of hardware chamber slot (zero-based)
Channelld	index of channel (zero-based)

Returns

example: return value -9 means that amps/volts per digit is in nano

gets the channel's unit name

Chamber⊷	index of hardware chamber slot (zero-based)
ld	
Channelld	index of channel (zero-based)

Returns

gets amps/volts per digit for specified chamber and channel

Parameters

Chamber⊷ Id	index of hardware chamber slot (zero-based)
Channelld	index of channel (zero-based)

Returns

```
11.121.3.23 GetUptimeSeconds() int32_t GetUptimeSeconds ( int32_t ChamberId )
```

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

```
11.121.3.24 GetVoltageClampControllerParam_D() uint32_t GetVoltageClampControllerParam_D ( int32_t ChamberId )
```

gets D value of PID controller;

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.25 GetVoltageClampControllerParam_I() uint32_t GetVoltageClampControllerParam_I (int32_t ChamberId)

gets I value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.26 GetVoltageClampControllerParam_P() uint32_t GetVoltageClampControllerParam_P (int32_t ChamberId)

gets P value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.27 IsChamberAvailable() bool IsChamberAvailable (int32_t *ChamberId*)

checks whether chamber amp is connected to slot

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.28 IsHighCurrentMode() bool IsHighCurrentMode (int32_t *ChamberId*)

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

11.121.3.29 IsInternalCalibrationFinished() bool IsInternalCalibrationFinished (int32_t ChamberId)

when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
Id	

Returns

returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED

Chamberld	index of hardware chamber slot (zero-based)	
ClampMode voltage or current clamp stimulation		

Returns

when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

sets clamp mode (voltage, current or open clamp)

Parameters

Chamberld	index of hardware chamber slot (zero-based)
NewClampMode	

```
11.121.3.32 SetEnablePulse() void SetEnablePulse (
    int32_t ChamberId,
    UssingClampModeEnumNet ClampMode,
    bool Enable )
```

enable pulse of given chamber and mode (voltage/current clamp) of this chamber

Parameters

Chamberld	index of hardware chamber slot (zero-based)	
ClampMode	voltage or current clamp stimulation	
Enable	when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level	

11.121.3.33 SetHighCurrentMode() void SetHighCurrentMode (int32_t *ChamberId*)

Parameters

Chamber⊷ Id

sets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

Chamberld	index of hardware chamber slot (zero-based)	
ClampMode	voltage or current clamp stimulation	
NewIdleOffset	unit: nA or mV	

Parameters

Chamberld	index of hardware chamber slot (zero-based)
NewLiquidResistance_Ohm	

```
11.121.3.36 SetLowCurrentMode() void SetLowCurrentMode (
    int32_t ChamberId )
```

Parameters

```
Chamber←
Id
```

```
11.121.3.37 SetPulse() void SetPulse (
    int32_t ChamberId,
    UssingClampModeEnumNet StgMode,
    int32_t NumberOfRepetitions,
    array< int >^ Amplitudes,
    array< int >^ Durations )
```

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

Chamberld	index of hardware chamber slot (zero-based); send pattern to connected amp	
StgMode	voltage or current clamp stimulation	
NumberOfRepetitions	number of repetitions for pulse pattern (-1 for infinite; n means pattern is applied n+1 times)	
Amplitudes	amplitude; unit in voltage clamp: mV; unit in current clamp: nA	
Durations	duration in 100us; CAUTION: first element is applied only one; auto-loop back to second element after last one	

11.121.3.38 SetVoltageClampControllerParam_D() void SetVoltageClampControllerParam_D (int32_t ChamberId, uint32_t D)

sets D value of PID controller;

Parameters

Chamber←		index of hardware chamber slot (zero-based)
ld		
D		useful range: 0700

11.121.3.39 SetVoltageClampControllerParam_I() void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)

sets I value of PID controller;

Parameters

Chamber⊷ Id	index of hardware chamber slot (zero-based)
1	useful range: 80000120000

11.121.3.40 SetVoltageClampControllerParam_P() void SetVoltageClampControllerParam_P (int32_t ChamberId, uint32_t P)

sets P value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	
Р	useful value: 130000

11.121.3.41 WaitForAllChambers() bool WaitForAllChambers ()

blocking call that waits for ALL chambers' boot-up calibration to complete

Returns

returns false when at least one chamber's calibration fails (e.g. timeout...)

blocking call that waits for chamber boot-up calibration to complete

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	


Returns

returns false when calibration fails (e.g. timeout...)

11.122 CWarnerValveControllerDeviceNet Class Reference

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

Inheritance diagram for CWarnerValveControllerDeviceNet:

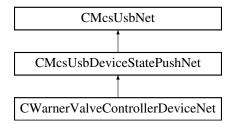

Public Member Functions

- delegate void OnGetValveActive (uint16 t valve, int valveActive)
- delegate void OnGetValveManualState (uint16_t valve, int32_t valveManualState)
- delegate void OnGetValveManualGroup (uint16_t valve, int32_t valveManualGroup)
- delegate void OnGetValveMode (uint16_t valve, WvcValveModeEnumNet ValveMode)
- delegate void OnGetAnalogThresholdLow (uint16_t valve, int32_t threshold)
- delegate void OnGetAnalogThresholdHigh (uint16_t valve, int32_t threshold)
- delegate void OnGetDigitalPortDirection (uint16_t port, PortDirectionEnumNet direction)

- delegate void OnlsValveDigitalInInverted (uint16_t valve, bool isInverted)
- delegate void OnGetValveDigitalInPort (uint16 t valve, uint32 t digitalInPort)
- delegate void OnlsDigitalOutPortInverted (uint16_t digitalOutPort, bool isInverted)
- delegate void OnGetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)
- delegate void OnlsValveOpen (uint16 t valve, bool valveOpen)
- delegate void OnlsValveOpenInDigitalMode (uint16_t valve, bool valveOpen)
- delegate void OnlsValveOpenInAnalogMode (uint16_t valve, bool valveOpen)
- delegate void OnGetAnalogVoltage (int32_t voltage)
- delegate void OnTableEntryChanged (uint16 t tableNumber)
- delegate void OnGetTableNamebyIndex (uint16_t tableNumber, String^ tableName)
- delegate void OnGetActiveRunningTableNumber (uint32 t tableNumber)
- delegate void OnGetCurrentNumberOfValves (int32_t numberOfValves)
- delegate void OnGetValveBoardRevision (uint32_t revision)
- delegate void OnGetValveLedOn (bool ledon)
- delegate void OnGetDisplayMode (WvcDisplayModeEnumNet DisplayMode)
- CWarnerValveControllerDeviceNet ()

Initializes a new instance of the CWarnerValveControllerDeviceNet class.

- virtual ~CWarnerValveControllerDeviceNet ()
- !CWarnerValveControllerDeviceNet ()
- int GetValveActive (uint16_t valve)

Gets the valve active/inactive state

void SetValveActive (uint16 t valve, int valveActive)

Sets the valve active/inactive state

uint32 t GetValvesActiveMap ()

Gets the valves active/inactive states

void SetValvesActiveMap (uint32_t valvesActive)

Sets the valve active/inactive state

• int32 t GetValveManualState (uint16 t valve)

Gets the valve manual on/off state

void SetValveManualState (uint16_t valve, int32_t valveManualState)

Sets the valve manual on/off state

uint32 t GetValvesManualStateMap ()

Gets the valves manual on/off states

void SetValvesManualStateMap (uint32_t valveaManualState)

Sets the valve manual on/off state

int32 t GetValveManualGroup (uint16 t valve)

Gets the valve manual group

void SetValveManualGroup (uint16_t valve, int32_t valveManualGroup)

Sets the valve manual group

WvcValveModeEnumNet GetValveMode (uint16 t valve)

Reads the valve mode

• void SetValveMode (uint16_t valve, WvcValveModeEnumNet ValveMode)

Writes the valve mode

int32_t GetAnalogThresholdLow (uint16_t valve)

Gets the lower threshold for the analog in port per valve

void SetAnalogThresholdLow (uint16_t valve, int32_t threshold)

Sets the lower threshold for the analog in port per valve

int32_t GetAnalogThresholdHigh (uint16_t valve)

Gets the upper threshold for the analog in port per valve

void SetAnalogThresholdHigh (uint16_t valve, int32_t threshold)

Sets the upper threshold for the analog in port per valve

PortDirectionEnumNet GetDigitalPortDirection (uint16_t port)

Gets the direction of a digital port

void SetDigitalPortDirection (uint16 t port, PortDirectionEnumNet direction)

Sets the direction of a digital port

bool IsValveDigitalInInverted (uint16_t valve)

Is digital in inverted

void SetValveDigitalInInvert (uint16_t valve, bool isInverted)

Invert digital in

uint32 t GetValveDigitalInPort (uint16 t valve)

Gets the number of the digital in port which is mapped to a valve

void SetValveDigitalInPort (uint16_t valve, uint32_t digitalInPort)

Map a digital in port to a valve

bool IsDigitalOutPortInverted (uint16_t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

void SetDigitalOutPortInvert (uint16_t digitalOutPort, bool isInverted)

Map a valve to a digital out port

uint32 t GetDigitalOutPortValve (uint16 t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

void SetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)

Map a valve to a digital out port

· void SetDefault ()

Sets the settings of the valve controller to default

bool IsValveOpen (uint16 t valve)

Is valve open

bool IsValveOpenInDigitalMode (uint16_t valve)

True, if the valve would be open when the device is in digital mode

bool IsValveOpenInAnalogMode (uint16_t valve)

True, if the valve would be open when the device is in analog mode

int32_t GetAnalogVoltage ()

Reads the voltage on the analog in port

void GetValveTableEntry (uint16_t valve, uint16_t index, [System::Runtime::InteropServices::Out]uint32_t% duration, [System::Runtime::InteropServices::Out]bool% state)

Read an entry from the valve protocol table

void SetValveTableEntry (uint16_t valve, uint16_t index, uint32_t duration, bool state)

Write an entry to the valve protocol table

void ClearValveTable (uint16_t valve)

Clear the valve protocol table

void LoadValveTable ()

Load the current table from permanent memory

void StoreValveTable ()

Store the current table in permanent memory

String \(^\) GetTableNamebyIndex (uint16_t tableNumber)

Get the name of a protocol table

String \(^\) GetTableName ()

Get the name of the current protocol table

void SetTableName (String[^] tableName)

Set the name of the current protocol table

uint32_t GetActiveRunningTableNumber ()

Gets the number of the table that is active for running

void SetActiveRunningTableNumber (uint32 t tableNumber)

Sets the number of the tanle that is active for running

uint32_t GetCurrentEditTableNumber ()

Gets the number of the table that is current for editing

void SetCurrentEditTableNumber (uint32 t tableNumber)

Sets the number of the table that is current for editing

void ClearTableName ()

Clear the name of current protocol table

void SetTableStep (uint16_t valve, int32_t steps)

Skips the table protocol for a valve by steps

void SetTableStepAll (int32 t steps)

Skips the table protocol for all valves by steps

• int32_t GetTotalNumberOfValves ()

Get the total number of valves in the system

• int32 t GetTotalNumberOfDigitalPorts ()

Get the total number of digital ports in the system

• int32 t GetTotalTableSize ()

Get the total table size in the system

int32 t GetTotalNumberOfTables ()

Get the total number of tables in the system

• int32 t GetCurrentNumberOfValves ()

Get the current number of valves connected to the system

uint32 t GetValveBoardRevision ()

Gets the revision code of the valve board

bool GetValveLedOn ()

Gets the LED state of the valve board

· void SetValveLedOn (bool ledon)

Gets the LED state of the valve board

WvcDisplayModeEnumNet GetDisplayMode ()

Reads the display mode

void SetDisplayMode (WvcDisplayModeEnumNet DisplayMode, int32_t lockTimeMs)

Writes the display mode

• String ^ GetValveBoardRevisionString ()

Gets the revision name of the valve board

Events

• OnGetValveActive^ GetValveActiveEvent [add, remove, raise]

Event fires when the valve state for the valve number has changed

• OnGetValveManualState^ GetValveManualStateEvent [add, remove, raise]

Event fires when the manual valve state for the valve number has changed

• OnGetValveManualGroup GetValveManualGroupEvent [add, remove, raise]

Event fires when the manual valve group for the valve number has changed

OnGetValveMode[^] GetValveModeEvent [add, remove, raise]

Event fires when the valve mode for the valve number has changed

• OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent [add, remove, raise]

Event fires when the threshold in mV for the valve number has changed

• OnGetAnalogThresholdHigh^ GetAnalogThresholdHighEvent [add, remove, raise]

Event fires when the threshold in mV for the valve number has changed

• OnGetDigitalPortDirection GetDigitalPortDirectionEvent [add, remove, raise]

Event fires when the direction for the port number has changed

• OnlsValveDigitalInInverted^ IsValveDigitalInInvertedEvent [add, remove, raise]

Event fires when is inverted for the valve number has changed

 $\bullet \ \ On Get Valve Digital In Port ^{\wedge} \ Get Valve Digital In Port Event \quad [\ \texttt{add, remove, raise}]$

Event fires when the digital in port for the valve number has changed

• OnlsDigitalOutPortInverted^ IsDigitalOutPortInvertedEvent [add, remove, raise]

Event fires when is inverted for the digital out port has changed

OnGetDigitalOutPortValve[^] GetDigitalOutPortValveEvent [add, remove, raise]

Event fires when the valve number for the digital out port has changed

• OnlsValveOpen^ lsValveOpenEvent [add, remove, raise]

Event fires when is open for the valve number has changed

• OnlsValveOpenInDigitalMode^ IsValveOpenInDigitalModeEvent [add, remove, raise]

Event fires when is open for the valve number has changed

• OnlsValveOpenInAnalogMode^ IsValveOpenInAnalogModeEvent [add, remove, raise]

Event fires when is open for the valve number has changed

OnGetAnalogVoltage[^] GetAnalogVoltageEvent [add, remove, raise]

Event fires when the voltage in mV has changed

• OnTableEntryChanged^ TableEntryChangedEvent [add, remove, raise]

Event fires when an entry of a table changed

• OnGetTableNamebyIndex^ GetTableNamebyIndexEvent [add, remove, raise]

Event fires when the name of the table for the table number has changed

• OnGetActiveRunningTableNumber^ GetActiveRunningTableNumberEvent [add, remove, raise] Event fires when the table number has changed

• OnGetCurrentNumberOfValves^ GetCurrentNumberOfValvesEvent [add, remove, raise]

Event fires when the number of valves has changed

• OnGetValveBoardRevision^ GetValveBoardRevisionEvent [add, remove, raise]

Event fires when the revision code has changed

• $OnGetValveLedOn^{\wedge} GetValveLedOnEvent$ [add, remove, raise]

Event fires when the LED state has changed

• OnGetDisplayMode^ GetDisplayModeEvent [add, remove, raise]

Event fires when the display mode has changed

Additional Inherited Members

11.122.1 Detailed Description

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

11.122.2 Constructor & Destructor Documentation

11.122.2.1 CWarnerValveControllerDeviceNet() CWarnerValveControllerDeviceNet ()

Initializes a new instance of the CWarnerValveControllerDeviceNet class.

11.122.2.2 ~CWarnerValveControllerDeviceNet() virtual ~CWarnerValveControllerDeviceNet () [virtual]

 $\textbf{11.122.2.3} \quad \textbf{"!CWarnerValveControllerDeviceNet()} \quad \textbf{!CWarnerValveControllerDeviceNet ()}$

11.122.3 Member Function Documentation

11.122.3.1 ClearTableName() void ClearTableName ()

Clear the name of current protocol table

```
11.122.3.2 ClearValveTable() void ClearValveTable ( uint16_t valve )
```

Clear the valve protocol table

Parameters

11.122.3.3 GetActiveRunningTableNumber() uint32_t GetActiveRunningTableNumber ()

Gets the number of the table that is active for running

Returns

The table number

```
11.122.3.4 GetAnalogThresholdHigh() int32_t GetAnalogThresholdHigh ( uint16_t valve)
```

Gets the upper threshold for the analog in port per valve

Parameters

valve The valve number

Returns

The threshold in mV

11.122.3.5 GetAnalogThresholdLow() int32_t GetAnalogThresholdLow (uint16_t valve)

Gets the lower threshold for the analog in port per valve

Parameters

valve The valve number

Returns

The threshold in mV

11.122.3.6 GetAnalogVoltage() int32_t GetAnalogVoltage ()

Reads the voltage on the analog in port

Returns

The voltage in mV

11.122.3.7 GetCurrentEditTableNumber() uint32_t GetCurrentEditTableNumber ()

Gets the number of the table that is current for editing

Returns

The table number

$\textbf{11.122.3.8} \quad \textbf{GetCurrentNumberOfValves()} \quad \texttt{int32_t GetCurrentNumberOfValves ()} \\$

Get the current number of valves connected to the system

Returns

The number of valves

```
11.122.3.9 GetDigitalOutPortValve() uint32_t GetDigitalOutPortValve ( uint16_t digitalOutPort )
```

Gets the number of the valve which is mapped to a digital out port

digitalOutPort	The digital out port
----------------	----------------------

Returns

The valve number

11.122.3.10 **GetDigitalPortDirection()** PortDirectionEnumNet GetDigitalPortDirection (uint16_t port)

Gets the direction of a digital port

Parameters

```
port The port number
```

Returns

the direction

11.122.3.11 GetDisplayMode() WvcDisplayModeEnumNet GetDisplayMode ()

Reads the display mode

Returns

The display mode

11.122.3.12 GetTableName() String $^{\land}$ GetTableName ()

Get the name of the current protocol table

Returns

The name of the table

11.122.3.13 GetTableNamebyIndex() String $^{\land}$ GetTableNamebyIndex (uint16_t tableNumber)

Get the name of a protocol table

tableNumber T	he table number
---------------	-----------------

Returns

The name of the table

11.122.3.14 GetTotalNumberOfDigitalPorts() int32_t GetTotalNumberOfDigitalPorts ()

Get the total number of digital ports in the system

Returns

The number of digital ports

11.122.3.15 GetTotalNumberOfTables() int32_t GetTotalNumberOfTables ()

Get the total number of tables in the system

Returns

The number of tables

11.122.3.16 GetTotalNumberOfValves() int32_t GetTotalNumberOfValves ()

Get the total number of valves in the system

Returns

The number of valves

11.122.3.17 GetTotalTableSize() int32_t GetTotalTableSize ()

Get the total table size in the system

Returns

The table size

```
11.122.3.18 GetValveActive() int GetValveActive ( uint16\_t \ valve )
```

Gets the valve active/inactive state

valve The valve number

Returns

The valve state

11.122.3.19 GetValveBoardRevision() uint32_t GetValveBoardRevision ()

Gets the revision code of the valve board

Returns

The revision code

11.122.3.20 GetValveBoardRevisionString() String ^ GetValveBoardRevisionString ()

Gets the revision name of the valve board

Returns

The revision name

11.122.3.21 GetValveDigitalInPort() uint32_t GetValveDigitalInPort (uint16_t *valve*)

Gets the number of the digital in port which is mapped to a valve

Parameters

valve	The valve number

Returns

The digital in port

11.122.3.22 GetValveLedOn() bool GetValveLedOn ()

Gets the LED state of the valve board

Returns

The LED state

$\textbf{11.122.3.23} \quad \textbf{GetValveManualGroup()} \quad \texttt{int32_t GetValveManualGroup} \quad \textbf{(}$

uint16_t valve)

Gets the valve manual group

Parameters

valve The valve no	umber
--------------------	-------

Returns

The manual valve group

11.122.3.24 GetValveManualState() int32_t GetValveManualState (uint16_t valve)

Gets the valve manual on/off state

Parameters

valve	The valve number

Returns

The manual valve state

11.122.3.25 **GetValveMode()** WvcValveModeEnumNet GetValveMode (uint16_t valve)

Reads the valve mode

Parameters

valve The valve number	
------------------------	--

Returns

The valve mode

11.122.3.26 GetValvesActiveMap() uint32_t GetValvesActiveMap ()

Gets the valves active/inactive states

Returns

The valves states

11.122.3.27 GetValvesManualStateMap() uint32_t GetValvesManualStateMap ()

Gets the valves manual on/off states

Returns

The manual valves states

11.122.3.28 GetValveTableEntry() void GetValveTableEntry (

```
uint16_t valve,
uint16_t index,
[System::Runtime::InteropServices::Out] uint32_t% duration,
[System::Runtime::InteropServices::Out] bool% state )
```

Read an entry from the valve protocol table

Parameters

valve	The valve number
index	The index in the table
duration	the duration in ms
state	the state

11.122.3.29 IsDigitalOutPortInverted() bool IsDigitalOutPortInverted (uint16_t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

Parameters

digitalOutPort	The digital out port

Returns

is inverted

```
11.122.3.30 IsValveDigitalInInverted() bool IsValveDigitalInInverted ( uint16_t valve)
```

Is digital in inverted

Parameters

```
valve The valve number
```

Returns

is inverted

```
11.122.3.31 IsValveOpen() bool IsValveOpen ( uint16_t valve)
```

_

Is valve open

Parameters

valve The valve number

Returns

is open

11.122.3.32 IsValveOpenInAnalogMode() bool IsValveOpenInAnalogMode (uint16_t *valve*)

True, if the valve would be open when the device is in analog mode

Parameters

```
valve The valve number
```

Returns

is open

```
11.122.3.33 IsValveOpenInDigitalMode() bool IsValveOpenInDigitalMode ( uint16_t valve)
```

True, if the valve would be open when the device is in digital mode

Pa	ra	m	ρi	ŀρ	re

Returns

is open

```
11.122.3.34 LoadValveTable() void LoadValveTable ( )
```

Load the current table from permanent memory

```
11.122.3.35 OnGetActiveRunningTableNumber() delegate void OnGetActiveRunningTableNumber ( uint32_t tableNumber)
```

```
11.122.3.36 OnGetAnalogThresholdHigh() delegate void OnGetAnalogThresholdHigh ( uint16_t valve, int32_t threshold)
```

```
11.122.3.37 OnGetAnalogThresholdLow() delegate void OnGetAnalogThresholdLow ( uint16_t valve, int32_t threshold)
```

```
11.122.3.39 OnGetCurrentNumberOfValves() delegate void OnGetCurrentNumberOfValves ( int32_t numberOfValves )
```

```
11.122.3.40 OnGetDigitalOutPortValve() delegate void OnGetDigitalOutPortValve ( uint16_t digitalOutPort, uint32_t valve )
```

```
11.122.3.41 OnGetDigitalPortDirection() delegate void OnGetDigitalPortDirection (
             uint16_t port,
             PortDirectionEnumNet direction )
11.122.3.42 OnGetDisplayMode() delegate void OnGetDisplayMode (
             WvcDisplayModeEnumNet DisplayMode )
11.122.3.43 OnGetTableNamebyIndex() delegate void OnGetTableNamebyIndex (
             uint16_t tableNumber,
             String^{\wedge} tableName)
11.122.3.44 OnGetValveActive() delegate void OnGetValveActive (
             uint16_t valve,
             int valveActive )
11.122.3.45 OnGetValveBoardRevision() delegate void OnGetValveBoardRevision (
             uint32_t revision )
11.122.3.46 OnGetValveDigitalInPort() delegate void OnGetValveDigitalInPort (
             uint16_t valve,
             uint32_t digitalInPort )
11.122.3.47 OnGetValveLedOn() delegate void OnGetValveLedOn (
             bool ledon )
11.122.3.48 OnGetValveManualGroup() delegate void OnGetValveManualGroup (
             uint16_t valve,
             int32_t valveManualGroup )
11.122.3.49 OnGetValveManualState() delegate void OnGetValveManualState (
             uint16_t valve,
             int32_t valveManualState )
```

```
11.122.3.50 OnGetValveMode() delegate void OnGetValveMode (
              uint16_t valve,
              WvcValveModeEnumNet ValveMode )
11.122.3.51 OnlsDigitalOutPortInverted() delegate void OnIsDigitalOutPortInverted (
              uint16_t digitalOutPort,
              bool isInverted )
\textbf{11.122.3.52} \quad \textbf{OnlsValveDigitalInInverted()} \quad \texttt{delegate void OnlsValveDigitalInInverted ()}
              uint16_t valve,
              bool is Inverted )
11.122.3.53 OnlsValveOpen() delegate void OnIsValveOpen (
              uint16_t valve,
              bool valveOpen )
11.122.3.54 OnlsValveOpenInAnalogMode() delegate void OnIsValveOpenInAnalogMode (
              uint16_t valve,
              bool valveOpen )
\textbf{11.122.3.55} \quad \textbf{OnlsValveOpenInDigitalMode()} \quad \texttt{delegate void OnlsValveOpenInDigitalMode ()} \\
              uint16_t valve,
              bool valveOpen )
11.122.3.56 OnTableEntryChanged() delegate void OnTableEntryChanged (
              uint16_t tableNumber )
11.122.3.57 SetActiveRunningTableNumber() void SetActiveRunningTableNumber (
              uint32_t tableNumber )
Sets the number of the tanle that is active for running
Parameters
```

tableNumber

The table number

11.122.3.58 SetAnalogThresholdHigh() void SetAnalogThresholdHigh (uint16_t valve, int32_t threshold)

Sets the upper threshold for the analog in port per valve

Parameters

valve	The valve number
threshold	The threshold in mV

Sets the lower threshold for the analog in port per valve

Parameters

valve	The valve number
threshold	The threshold in mV

11.122.3.60 SetCurrentEditTableNumber() void SetCurrentEditTableNumber (uint32_t tableNumber)

Sets the number of the table that is current for editing

Parameters

11.122.3.61 SetDefault() void SetDefault ()

Sets the settings of the valve controller to default

Map a valve to a digital out port

digitalOutPort	The digital out port
isInverted	True if digital out is to be inverted

11.122.3.63 SetDigitalOutPortValve() void SetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)

Map a valve to a digital out port

Parameters

digitalOutPort	The digital out port
valve	The valve number

Sets the direction of a digital port

Parameters

port	The port number
direction	the direction

Writes the display mode

DisplayMode	The display mode
lockTimeMs	Locks the display for ms

```
11.122.3.66 SetTableName() void SetTableName ( String^{ } tableName )
```

Set the name of the current protocol table

Parameters

tableName	The name of the table
-----------	-----------------------

Skips the table protocol for a valve by steps

Parameters

valve	The valve number
steps	Number of steps

Skips the table protocol for all valves by steps

Parameters

steps	Number of steps

Sets the valve active/inactive state

valve	The valve number
valveActive	The valve state

```
11.122.3.70 SetValveDigitalInInvert() void SetValveDigitalInInvert ( uint16_t valve, bool isInverted )
```

Invert digital in

Parameters

valve	The valve number
isInverted	True if digital in is to be inverted

Map a digital in port to a valve

Parameters

valve	The valve number
digitalInPort	The digital in port

Gets the LED state of the valve board

Parameters

ledon The LED state

Sets the valve manual group

valve	The valve number
valveManualGroup	The manual valve group

Sets the valve manual on/off state

Parameters

valve	The valve number
valveManualState	The manual valve state

Writes the valve mode

Parameters

valve	The valve number
ValveMode	The valve mode

11.122.3.76 SetValvesActiveMap() void SetValvesActiveMap (uint32_t valvesActive)

Sets the valve active/inactive state

Parameters

valvesActive	The valves states
vaiveshelive	The valves states

11.122.3.77 SetValvesManualStateMap() void SetValvesManualStateMap (uint32_t valveaManualState)

Sets the valve manual on/off state

Parameters

```
valveaManualState The manual valves states
```


Write an entry to the valve protocol table

Parameters

valve	The valve number
index	The index in the table
duration	the duration in ms
state	the state

11.122.3.79 StoreValveTable() void StoreValveTable ()

Store the current table in permanent memory

11.122.4 Event Documentation

11.122.4.1 GetActiveRunningTableNumberEvent OnGetActiveRunningTableNumber^ GetActiveRunning← TableNumberEvent [add], [remove], [raise]

Event fires when the table number has changed

11.122.4.2 GetAnalogThresholdHighEvent OnGetAnalogThresholdHigh^ GetAnalogThresholdHighEvent [add], [remove], [raise]

Event fires when the threshold in mV for the valve number has changed

11.122.4.3 GetAnalogThresholdLowEvent OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent [add], [remove], [raise]

Event fires when the threshold in $\ensuremath{\mathsf{mV}}$ for the valve number has changed

11.122.4.4 GetAnalogVoltageEvent OnGetAnalogVoltage^ GetAnalogVoltageEvent [add], [remove], [raise]

Event fires when the voltage in mV has changed

11.122.4.5 GetCurrentNumberOfValvesEvent OnGetCurrentNumberOfValves^ GetCurrentNumberOf← ValvesEvent [add], [remove], [raise]

Event fires when the number of valves has changed

11.122.4.6 GetDigitalOutPortValveEvent OnGetDigitalOutPortValve^ GetDigitalOutPortValveEvent [add], [remove], [raise]

Event fires when the valve number for the digital out port has changed

11.122.4.7 GetDigitalPortDirectionEvent OnGetDigitalPortDirection^ GetDigitalPortDirectionEvent [add], [remove], [raise]

Event fires when the direction for the port number has changed

11.122.4.8 GetDisplayModeEvent OnGetDisplayMode^ GetDisplayModeEvent [add], [remove], [raise]

Event fires when the display mode has changed

 $\textbf{11.122.4.9} \quad \textbf{GetTableNamebyIndexEvent} \quad \texttt{OnGetTableNamebyIndex}^{\wedge} \quad \texttt{GetTableNamebyIndexEvent} \quad \texttt{[add],} \\ \texttt{[remove], [raise]}$

Event fires when the name of the table for the table number has changed

11.122.4.10 GetValveActiveEvent OnGetValveActive^ GetValveActiveEvent [add], [remove], [raise]

Event fires when the valve state for the valve number has changed

11.122.4.11 GetValveBoardRevisionEvent OnGetValveBoardRevision^ GetValveBoardRevisionEvent [add], [remove], [raise]

Event fires when the revision code has changed

11.122.4.12 GetValveDigitalInPortEvent OnGetValveDigitalInPort^ GetValveDigitalInPortEvent [add], [remove], [raise]

Event fires when the digital in port for the valve number has changed

11.122.4.13 GetValveLedOnEvent OnGetValveLedOn^ GetValveLedOnEvent [add], [remove], [raise]

Event fires when the LED state has changed

11.122.4.14 GetValveManualGroupEvent OnGetValveManualGroup^ GetValveManualGroupEvent [add], [remove], [raise]

Event fires when the manual valve group for the valve number has changed

11.122.4.15 GetValveManualStateEvent OnGetValveManualState^ GetValveManualStateEvent [add], [remove], [raise]

Event fires when the manual valve state for the valve number has changed

11.122.4.16 GetValveModeEvent OnGetValveMode^ GetValveModeEvent [add], [remove], [raise]

Event fires when the valve mode for the valve number has changed

11.122.4.17 IsDigitalOutPortInvertedEvent OnIsDigitalOutPortInverted[∧] IsDigitalOutPortInverted← Event [add], [remove], [raise]

Event fires when is inverted for the digital out port has changed

11.122.4.18 IsValveDigitalInInvertedEvent OnIsValveDigitalInInverted^ IsValveDigitalInInverted← Event [add], [remove], [raise]

Event fires when is inverted for the valve number has changed

11.122.4.19 IsValveOpenEvent OnIsValveOpen^ IsValveOpenEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.122.4.20 IsValveOpenInAnalogModeEvent OnIsValveOpenInAnalogMode^ IsValveOpenInAnalogMode← Event [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.122.4.21 IsValveOpenInDigitalModeEvent OnIsValveOpenInDigitalMode^ IsValveOpenInDigital← ModeEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed


11.122.4.22 TableEntryChangedEvent OnTableEntryChanged^ TableEntryChangedEvent [add], [remove], [raise]

Event fires when an entry of a table changed

11.123 CWarnerValveControllerDeviceTesterFunctionNet Class Reference

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

Inheritance diagram for CWarnerValveControllerDeviceTesterFunctionNet:

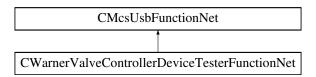

Public Member Functions

CWarnerValveControllerDeviceTesterFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pWarnerValveControllerDeviceTesterFunctionPointerContainer)

Initializes a new instance of the CWarnerValveControllerDeviceTesterFunctionNet class.

- CWarnerValveControllerDeviceTesterFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CWarnerValveControllerDeviceTesterFunctionNet ()
- !CWarnerValveControllerDeviceTesterFunctionNet ()
- void SetADC (uint32_t onoff)

Sets the ADC port of the tester

• uint32_t GetSync ()

Gets the output from the sync port

void SetTrigger (uint32 t trigger)

Sets the input to the trigger port

void SetTriggerSyncDirection (uint32_t direction)

Sets the direction of the trigger/sync test port

· uint32_t GetIO ()

Gets the output from the io ports

void SetIO (uint32_t io)

Sets the input to the io ports

· void SetIODirection (int32_t direction)

Sets the direction of the IO test ports

Additional Inherited Members

11.123.1 Detailed Description

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

11.123.2 Constructor & Destructor Documentation

Initializes a new instance of the CWarnerValveControllerDeviceTesterFunctionNet class.

```
11.123.2.2 CWarnerValveControllerDeviceTesterFunctionNet() [2/2] CWarnerValveControllerDeviceTesterFunctionNet (

CMcsUsbNet^ mcsusb )
```

```
11.123.2.3 ~CWarnerValveControllerDeviceTesterFunctionNet() virtual ~CWarnerValveControllerDeviceTesterFunctio
( ) [virtual]
```

```
11.123.2.4 "!CWarnerValveControllerDeviceTesterFunctionNet() !CWarnerValveControllerDeviceTesterFunctionNet
```

11.123.3 Member Function Documentation

```
11.123.3.1 GetIO() uint32_t GetIO ()
```

Gets the output from the io ports

Returns

The manual valves states

```
11.123.3.2 GetSync() uint32_t GetSync ()
```

Gets the output from the sync port

Returns

The sync state

```
11.123.3.3 SetADC() void SetADC ( uint32_t onoff )
```

Sets the ADC port of the tester

Parameters

```
onoff The port state
```

```
11.123.3.4 SetIO() void SetIO ( uint32_t io )
```

Sets the input to the io ports

Parameters

io The manual valves states

```
11.123.3.5 SetIODirection() void SetIODirection ( int32_t direction )
```

Sets the direction of the IO test ports

Parameters

```
direction The 16bit direction map: 1=IN 0=OUT
```

Sets the input to the trigger port

trigger	The trigger state
---------	-------------------

11.123.3.7 SetTriggerSyncDirection() void SetTriggerSyncDirection (uint32_t direction)


Sets the direction of the trigger/sync test port

Parameters

direction The direction: 1=IN 0=OUT

11.124 CWClassicFunctionNet Class Reference

Inheritance diagram for CWClassicFunctionNet:

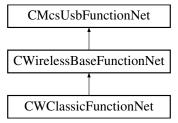

Public Member Functions

- CWClassicFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] wClassicFuntion
 — PointerContainer)
- CWClassicFunctionNet (CMcsUsbNet[^] mcsusb)
- uint32 t ResetChannelmap (unsigned int virtualDevice)
- uint32 t SetChannelmap (unsigned char position, unsigned char channel, unsigned int Device)
- void SetHWSelectedChannels (array< bool >^ channels, unsigned int Device)
- void SetRFLostBehaviour (uint8 t stoponfailure, unsigned int Device)
- void SetHeadstageOnOff (uint16_t onoff)
- USHORT GetHeadstageOnOff ()
- void SetRFFrequencyHeadstage (uint8_t receiver_nb, unsigned short frequency)
- unsigned short GetRFFrequencyHeadstage (uint8_t receiver_nb)
- · void SetRFFrequencyReceiver (uint8 t receiver nb, uint8 t configuration, unsigned short frequency)
- void SetRFFrequencyReceiverEeprom (uint8_t receiver_nb, uint8_t configuration, unsigned short frequency)
- unsigned short GetRFFrequencyReceiver (uint8 t receiver nb, uint8 t configuration)
- void SetSerialNumberHeadstage (unsigned short number)
- unsigned short GetSerialNumberHeadstage ()
- void SetSelectedHeadstage (uint8_t number)
- uint8_t GetSelectedHeadstage ()
- void ScanForHeadstages ()
- uint8_t GetScanHeadstagesResult (int max_wait_for_ms)

- void SetFilterParametersHeadstage (unsigned short index, array< int >^ buffer)
- array< int > ^ GetFilterParametersHeadstage (unsigned short index)
- bool GetHasRedLedHeadstage ()
- · void SetHasChecksum (unsigned int has, unsigned int Device)
- unsigned int GetHasChecksum (unsigned int Device)
- void SetResetFilter (unsigned int reset, unsigned int Device)
- unsigned int GetResetFilter (unsigned int Device)
- void SetWPAType (unsigned short type, unsigned int Device)
- · unsigned short GetWPAType (unsigned int Device)
- void SetWPADebugMode (unsigned int mode, unsigned int Device)
- unsigned int GetWPADebugMode (unsigned int Device)
- void SetRFPower (unsigned short power)
- unsigned short GetRFPower ()
- unsigned int GetRFConnectionStatus ()

Additional Inherited Members

11.124.1 Constructor & Destructor Documentation

11.124.2 Member Function Documentation

11.124.2.3 GetHasRedLedHeadstage() bool GetHasRedLedHeadstage ()


```
11.124.2.4 GetHeadstageOnOff() USHORT GetHeadstageOnOff ()
11.124.2.5 GetResetFilter() unsigned int GetResetFilter (
             unsigned int Device )
\textbf{11.124.2.6} \quad \textbf{GetRFConnectionStatus()} \quad \textbf{unsigned int GetRFConnectionStatus ()}
11.124.2.7 GetRFFrequencyHeadstage() unsigned short GetRFFrequencyHeadstage (
             uint8_t receiver_nb )
11.124.2.8 GetRFFrequencyReceiver() unsigned short GetRFFrequencyReceiver (
             uint8_t receiver_nb,
             uint8_t configuration )
11.124.2.9 GetRFPower() unsigned short GetRFPower ()
11.124.2.10 GetScanHeadstagesResult() uint8_t GetScanHeadstagesResult (
             int max\_wait\_for\_ms )
11.124.2.11 GetSelectedHeadstage() uint8_t GetSelectedHeadstage ( )
11.124.2.12 GetSerialNumberHeadstage() unsigned short GetSerialNumberHeadstage ( )
11.124.2.13 GetWPADebugMode() unsigned int GetWPADebugMode (
             unsigned int Device )
```

```
\textbf{11.124.2.14} \quad \textbf{GetWPAType()} \quad \texttt{unsigned short GetWPAType (}
             unsigned int Device )
11.124.2.15 ResetChannelmap() uint32_t ResetChannelmap (
             unsigned int virtualDevice )
11.124.2.16 ScanForHeadstages() void ScanForHeadstages ()
11.124.2.17 SetChannelmap() uint32_t SetChannelmap (
             unsigned char position,
             unsigned char channel,
             unsigned int Device )
11.124.2.18 SetFilterParametersHeadstage() void SetFilterParametersHeadstage (
             unsigned short index,
             array< int >^{\wedge} buffer )
11.124.2.19 SetHasChecksum() void SetHasChecksum (
             unsigned int has,
             unsigned int Device )
11.124.2.20 SetHeadstageOnOff() void SetHeadstageOnOff (
             uint16_t onoff )
11.124.2.21 SetHWSelectedChannels() void SetHWSelectedChannels (
             array< bool >^{\wedge} channels,
             unsigned int Device )
11.124.2.22 SetResetFilter() void SetResetFilter (
             unsigned int reset,
             unsigned int Device )
```

```
\textbf{11.124.2.23} \quad \textbf{SetRFF} \textbf{requencyHeadstage()} \quad \textbf{void SetRFF} \textbf{requencyHeadstage ()}
              uint8_t receiver_nb,
              unsigned short frequency )
11.124.2.24 SetRFFrequencyReceiver() void SetRFFrequencyReceiver (
              uint8_t receiver_nb,
              uint8_t configuration,
              unsigned short frequency )
11.124.2.25 SetRFFrequencyReceiverEeprom() void SetRFFrequencyReceiverEeprom (
              uint8_t receiver_nb,
              uint8_t configuration,
              unsigned short frequency )
11.124.2.26 SetRFLostBehaviour() void SetRFLostBehaviour (
              uint8_t stoponfailure,
              unsigned int Device )
11.124.2.27 SetRFPower() void SetRFPower (
              unsigned short power )
11.124.2.28 SetSelectedHeadstage() void SetSelectedHeadstage (
              uint8_t number )
11.124.2.29 SetSerialNumberHeadstage() void SetSerialNumberHeadstage (
              unsigned short number )
{\bf 11.124.2.30} \quad {\bf SetWPADebugMode()} \quad {\tt void \ SetWPADebugMode\ ()}
              unsigned int mode,
              unsigned int Device )
11.124.2.31 SetWPAType() void SetWPAType (
              unsigned short type,
              unsigned int Device )
```

11.125 CWirelessBaseFunctionNet Class Reference

Inheritance diagram for CWirelessBaseFunctionNet:

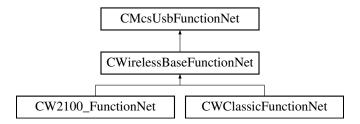

Public Member Functions

• CWirelessBaseFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ mcsusbfunction)

Static Public Member Functions

• static String \(^\) CreateWirelessHeadstageSerialNumberString (unsigned short ID)

Additional Inherited Members

11.125.1 Constructor & Destructor Documentation

```
11.125.1.1 CWirelessBaseFunctionNet() CWirelessBaseFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ mcsusbfunction)
```

11.125.2 Member Function Documentation

11.126 DeviceIdNet Struct Reference

Device Id.

Public Member Functions

- DeviceIdNet ()
- DeviceIdNet (VendorIdEnumNet vendor, ProductIdEnumNet product, int bcd, McsBusTypeEnumNet bustype)
- DeviceIdNet (DeviceIdNet% deviceId)
- DeviceIdNet operator= (DeviceIdNet% deviceId)

Public Attributes

- VendorldEnumNet IdVendor
- ProductIdEnumNet IdProduct
- int BcdDevice
- McsBusTypeEnumNet BusType

11.126.1 Detailed Description

Device Id.

11.126.2 Constructor & Destructor Documentation

```
11.126.2.1 DeviceIdNet() [1/3] DeviceIdNet ( )
```

```
11.126.2.2 DeviceIdNet() [2/3] DeviceIdNet (

VendorIdEnumNet vendor,
```

ProductIdEnumNet product, int bcd,

 ${\tt McsBusTypeEnumNet}\ bustype\)$

```
11.126.2.3 DeviceIdNet() [3/3] DeviceIdNet (

DeviceIdNet% deviceId)
```

11.126.3 Member Function Documentation

```
11.126.3.1 operator=() DeviceIdNet operator= (
DeviceIdNet% deviceId)
```

11.126.4 Member Data Documentation

```
11.126.4.1 BcdDevice int BcdDevice
```

```
11.126.4.2 BusType McsBusTypeEnumNet BusType
```

```
11.126.4.3 IdProduct ProductIdEnumNet IdProduct
```

```
11.126.4.4 IdVendor VendorIdEnumNet IdVendor
```

11.127 DigitalSource< digitalsourceenum > Class Template Reference

Public Member Functions

- DigitalSource (digitalsourceenum source)
- int MaxBitNumber ()

Static Public Member Functions

• static int size ()

Properties

• digitalsourceenum Source [get]

11.127.1 Constructor & Destructor Documentation

```
11.127.1.1 DigitalSource() DigitalSource (
digitalsourceenum source)
```

11.127.2 Member Function Documentation

```
11.127.2.1 MaxBitNumber() int MaxBitNumber ( )
```

```
11.127.2.2 size() static int size ( ) [static]
```

11.127.3 Property Documentation

11.127.3.1 Source digitalsourceenum Source [get]

11.128 DriverVersionNet Class Reference

Class gives firmware versions of the device's firmware destinations.

Public Member Functions

DriverVersionNet ()

Contructor.

∼DriverVersionNet ()

Destructor.

unsigned int GetStatus (CFirmwareDestinationNet dest)

Get status of firmware destination.

• unsigned int GetStatus (unsigned int index)

Get status of firmware destination.

• unsigned int GetVersionInt (CFirmwareDestinationNet dest)

Get the version number of firmware destination (major in high word, minor in low word)

unsigned int GetVersionInt (unsigned int index)

Get the version number of firmware destination (major in high word, minor in low word)

unsigned int GetMajor (CFirmwareDestinationNet dest)

Get the major version number of firmware destination.

unsigned int GetMajor (unsigned int index)

Get the major version number of firmware destination.

unsigned int GetMinor (CFirmwareDestinationNet dest)

Get the minor version number of firmware destination.

unsigned int GetMinor (unsigned int index)

Get the minor version number of firmware destination.

• unsigned int GetNumEntries ()

Get the number of available firmware destinations.

String \(^\) GetVersionString (CFirmwareDestinationNet dest)

Get the version as a string in the format Major. Minor.

String \(^\) GetVersionString (unsigned int index)

Get the version as a string in the format Major. Minor.

CFirmwareDestinationNet GetDestinationCode (unsigned int index)

Get CFirmwareDestinationNet.

String \(^\) GetDestinationName (CFirmwareDestinationNet dest)

Get firmware destination name.

String \(^\) GetDestinationName (unsigned int index)

Get firmware destination name.

Static Public Member Functions

static String \(^\) DriverVersionNet::FormatVersion (unsigned int v)

11.128.1 Detailed Description

Class gives firmware versions of the device's firmware destinations.

11.128.2 Constructor & Destructor Documentation

```
11.128.2.1 DriverVersionNet() DriverVersionNet ( )
```

Contructor.

```
11.128.2.2 ~DriverVersionNet() ~DriverVersionNet ()
```

Destructor.

11.128.3 Member Function Documentation

```
11.128.3.1 DriverVersionNet::FormatVersion() static String ^{\land} DriverVersionNet::FormatVersion ( unsigned int v ) [static]
```

```
11.128.3.2 GetDestinationCode() CFirmwareDestinationNet GetDestinationCode ( unsigned int index )
```

Get CFirmwareDestinationNet.

```
index by index of firmware destination
```

```
11.128.3.3 GetDestinationName() [1/2] String ^{\wedge} GetDestinationName ( CFirmwareDestinationNet dest )
```

Get firmware destination name.

dest by CFirmwareDestionationNet

11.128.3.4 GetDestinationName() [2/2] String ^ GetDestinationName (unsigned int index)

Get firmware destination name.

Parameters

index by index of firmware destination

11.128.3.5 GetMajor() [1/2] unsigned int GetMajor (CFirmwareDestinationNet dest)

Get the major version number of firmware destination.

Parameters

dest	by CFirmwareDestionationNet
------	-----------------------------

```
11.128.3.6 GetMajor() [2/2] unsigned int GetMajor ( unsigned int index )
```

Get the major version number of firmware destination.

Parameters

index by index of firmware destination

```
11.128.3.7 GetMinor() [1/2] unsigned int GetMinor (
CFirmwareDestinationNet dest)
```

Get the minor version number of firmware destination.

dest	by CFirmwareDestionationNet

```
11.128.3.8 GetMinor() [2/2] unsigned int GetMinor ( unsigned int index )
```

Get the minor version number of firmware destination.

Parameters

index	by index of firmware destination
-------	----------------------------------

11.128.3.9 GetNumEntries() unsigned int GetNumEntries ()

Get the number of available firmware destinations.

```
11.128.3.10 GetStatus() [1/2] unsigned int GetStatus (
CFirmwareDestinationNet dest)
```

Get status of firmware destination.

Parameters

dest by CFirmwareDestionationNet

```
11.128.3.11 GetStatus() [2/2] unsigned int GetStatus ( unsigned int index )
```

Get status of firmware destination.

Parameters

index by index of firmware destination

```
11.128.3.12 GetVersionInt() [1/2] unsigned int GetVersionInt (
CFirmwareDestinationNet dest)
```

Get the version number of firmware destination (major in high word, minor in low word)

dest by CFirmwareDestionationNet

11.128.3.13 GetVersionInt() [2/2] unsigned int GetVersionInt (unsigned int index)

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

index by index of firmware destination

```
11.128.3.14 GetVersionString() [1/2] String ^{\land} GetVersionString ( CFirmwareDestinationNet dest )
```

Get the version as a string in the format Major.Minor.

Parameters

dest by CFirmwareDestionationNet

```
11.128.3.15 GetVersionString() [2/2] String ^{\land} GetVersionString ( unsigned int index )
```

Get the version as a string in the format Major.Minor.

Parameters

index by index of firmware

11.129 Firmware Destination Names Class Reference

Static Public Attributes

```
    static String \(^\text{DSP} = \text{gcnew String( "DSP" )}\)
    static String \(^\text{USB} = \text{gcnew String( "USB" )}\)
```

• static String $^{\wedge}$ MCU1 = gcnew String("MCU1")

- static String \(^\text{ Bootstrap} = \text{gcnew String("Bootstrap")}\)
- static String \(^{\text{MCSBUS1}} = \text{gcnew String("McsBus1")}\)

```
    static String \(^{\text{MCSBUS2}} = \text{gcnew String( "McsBus2" )}\)

    static String \(^{\text{MCSBUS3}} = \text{gcnew String( "McsBus3" )}\)

    static String \(^\text{MCSBUS4} = \text{gcnew String( "McsBus4" )}\)

    static String \(^\text{MCSBUS5} = \text{gcnew String( "McsBus5" )}\)

    static String \(^{\text{MCSBUS6}} = \text{gcnew String( "McsBus6" )}\)

    static String \(^{\text{MCSBUS7}} = \text{gcnew String( "McsBus7" )}\)

    static String \(^{\text{MCSBUS8}} = \text{gcnew String( "McsBus8" )}\)

    static String \(^{\text{MCSBUS9}} = \text{gcnew String( "McsBus9")}\)

    static String \(^{\text{MCSBUS10}} = \text{gcnew String( "McsBus10" )}\)

    static String \(^\text{MCSBUS11} = \text{gcnew String("McsBus11")}\)

    static String ^ MCSBUS12 = gcnew String( "McsBus12" )

    static String \(^{\text{MCSBUS13}}\) = gcnew String( "McsBus13" )

    static String \(^\text{BUS1_MCSBUS1}\) = gcnew String("Bus1McsBus1")

    static String \(^\text{BUS1_MCSBUS2} = \text{gcnew String("Bus1McsBus2")}\)

    static String \(^\text{PIC} = \text{gcnew String("PIC")}\)

    static String \(^\text{PIC2} = \text{gcnew String("PIC2")}\)

    static String \(^\text{PIC3} = \text{gcnew String("PIC3")}\)

    static String \(^\text{PIC4} = \text{gcnew String("PIC4")}\)

• static String ^ Altera = gcnew String( "Altera" )

    static String \(^\text{FPGA2} = \text{gcnew String( "FPGA2" )}\)

    static String \(^\text{FPGA3} = \text{gcnew String("FPGA3")}\)

    static String \(^{\text{FPGA4}} = \text{gcnew String("FPGA4")}\)

    static String \(^{\text{FPGA5}} = \text{gcnew String( "FPGA5" )}\)

    static String \(^{\text{FPGA6}} = \text{gcnew String( "FPGA6" )}\)
```

11.129.1 Member Data Documentation

```
11.129.1.1 Altera String ^ Altera = gcnew String( "Altera" ) [static]

11.129.1.2 Bootstrap String ^ Bootstrap = gcnew String( "Bootstrap" ) [static]

11.129.1.3 BUS1_MCSBUS1 String ^ BUS1_MCSBUS1 = gcnew String( "Bus1McsBus1" ) [static]

11.129.1.4 BUS1_MCSBUS2 String ^ BUS1_MCSBUS2 = gcnew String( "Bus1McsBus2" ) [static]

11.129.1.5 DSP String ^ DSP = gcnew String( "DSP" ) [static]
```

```
11.129.1.6 FPGA2 String ^ FPGA2 = gcnew String( "FPGA2" ) [static]
11.129.1.7 FPGA3 String ^{\wedge} FPGA3 = gcnew String( "FPGA3" ) [static]
11.129.1.8 FPGA4 String ^{\wedge} FPGA4 = gcnew String( "FPGA4" ) [static]
11.129.1.9 FPGA5 String ^{\wedge} FPGA5 = gcnew String( "FPGA5" ) [static]
11.129.1.10 FPGA6 String ^ FPGA6 = gcnew String( "FPGA6" ) [static]
11.129.1.11 MCSBUS1 String ^ MCSBUS1 = gcnew String( "McsBus1" ) [static]
11.129.1.12 MCSBUS10 String ^ MCSBUS10 = gcnew String( "McsBus10") [static]
11.129.1.13 MCSBUS11 String ^ MCSBUS11 = gcnew String( "McsBus11" ) [static]
\textbf{11.129.1.14} \quad \textbf{MCSBUS12} \quad \texttt{String} \; \wedge \; \texttt{MCSBUS12} \; = \; \texttt{gcnew} \; \; \texttt{String("McsBus12")} \; \; [\; \texttt{static}]
11.129.1.15 MCSBUS13 String ^ MCSBUS13 = gcnew String( "McsBus13" ) [static]
11.129.1.16 MCSBUS2 String ^{\land} MCSBUS2 = gcnew String( "McsBus2" ) [static]
```


```
11.129.1.17 MCSBUS3 String ^{\wedge} MCSBUS3 = gcnew String( "McsBus3" ) [static]
11.129.1.18 MCSBUS4 String ^{\land} MCSBUS4 = gcnew String( "McsBus4" ) [static]
11.129.1.19 MCSBUS5 String ^{\wedge} MCSBUS5 = gcnew String( "McsBus5" ) [static]
11.129.1.20 MCSBUS6 String ^ MCSBUS6 = gcnew String( "McsBus6" ) [static]
11.129.1.21 MCSBUS7 String ^ MCSBUS7 = gcnew String( "McsBus7" ) [static]
11.129.1.22 MCSBUS8 String ^ MCSBUS8 = gcnew String( "McsBus8") [static]
11.129.1.23 MCSBUS9 String ^ MCSBUS9 = gcnew String( "McsBus9" ) [static]
11.129.1.24 MCU1 String ^ MCU1 = gcnew String( "MCU1" ) [static]
11.129.1.25 PIC String ^{\land} PIC = gcnew String( "PIC" ) [static]
11.129.1.26 PIC2 String ^{\wedge} PIC2 = gcnew String( "PIC2" ) [static]
11.129.1.27 PIC3 String ^{\wedge} PIC3 = gcnew String( "PIC3") [static]
```

```
11.129.1.28 PIC4 String ^ PIC4 = gcnew String( "PIC4" ) [static]
```

11.129.1.29 USB String ^ USB = gcnew String("USB") [static]

11.130 HeadStageIDType Class Reference

Inheritance diagram for HeadStageIDType:

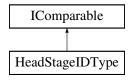

Public Types

enum HeadstageTypeEnum {
 Unknown,
 MeasuringOnly,
 OpticalStimulation,
 ElectricalStimulation }

Public Member Functions

- HeadStageIDType (unsigned int entry, CW2100_FunctionNet[^] device)
- virtual System::String ^ ToString () override
- virtual bool Equals (Object[^] obj) override
- virtual Int32 CompareTo (Object[^] obj)

Properties

- bool Valid [get]
- unsigned int Entry [get]
- unsigned short ID [get]
- System::String^ SN [get]
- unsigned int TypeValue [get]
- System::String^ Type [get]
- HeadstageTypeEnum HeadstageType [get]
- System::String^ UserDefinedName [get]
- int NumberOfAnalogChannels [get]
- int NumberOfStimulationChannels [get]
- W2100_StimulusParametersNet^{\(\Lambda\)} StimulusParameters [get]
- bool HasIMU [get]
- bool W16lsW14 [get]
- bool HasOptoCurrentMessurement [get]

11.130.1 Member Enumeration Documentation

11.130.1.1 HeadstageTypeEnum enum HeadstageTypeEnum [strong]

Enumerator

Unknown	
MeasuringOnly	
OpticalStimulation	
ElectricalStimulation	

11.130.2 Constructor & Destructor Documentation

11.130.3 Member Function Documentation

```
11.130.3.1 CompareTo() virtual Int32 CompareTo (
Object^ obj ) [virtual]
```

```
11.130.3.2 Equals() virtual bool Equals (
Object^ obj ) [override], [virtual]
```

```
11.130.3.3 ToString() virtual System::String ^ ToString ( ) [override], [virtual]
```

11.130.4 Property Documentation

11.130.4.1 Entry unsigned int Entry [get]

11.130.4.2 HasIMU bool HasIMU [get]

```
11.130.4.3 HasOptoCurrentMessurement bool HasOptoCurrentMessurement [get]
11.130.4.4 HeadstageType HeadstageTypeEnum HeadstageType [get]
11.130.4.5 ID unsigned short ID [get]
\textbf{11.130.4.6} \quad \textbf{NumberOfAnalogChannels} \quad \texttt{int NumberOfAnalogChannels} \quad \texttt{[get]}
11.130.4.7 NumberOfStimulationChannels int NumberOfStimulationChannels [get]
11.130.4.8 SN System:: String^ SN [get]
11.130.4.9 StimulusParameters W2100_StimulusParametersNet^ StimulusParameters [get]
11.130.4.10 Type System:: String^ Type [get]
11.130.4.11 TypeValue unsigned int TypeValue [get]
11.130.4.12 UserDefinedName System:: String^ UserDefinedName [get]
11.130.4.13 Valid bool Valid [get]
11.130.4.14 W16IsW14 bool W16IsW14 [get]
```

11.131 HeadstageIDTypeObject Class Reference

Public Member Functions

- HeadstageIDTypeObject (HeadStageIDType^ idType)
- virtual String \(^{\text{ToString}}\) () override
- virtual bool Equals (Object[∧] obj) override
- · virtual int GetHashCode () override

Public Attributes

- HeadStageIDType ^ _IdType
- String ^ _AdditionalText

Properties

- String^ AdditionalText [get, set]

11.131.1 Constructor & Destructor Documentation

```
11.131.1.1 HeadstageIDTypeObject() HeadstageIDTypeObject ( HeadStageIDType^ idType )
```

11.131.2 Member Function Documentation

```
11.131.2.1 Equals() virtual bool Equals (
          Object^ obj ) [override], [virtual]
```

```
11.131.2.2 GetHashCode() virtual int GetHashCode ( ) [override], [virtual]
```

```
11.131.2.3 ToString() virtual String ^{\wedge} ToString ( ) [override], [virtual]
```

11.131.3 Member Data Documentation

```
\textbf{11.131.3.1} \quad \textbf{\_AdditionalText} \quad \texttt{String} \; \land \; \textbf{\_AdditionalText}
```

```
11.131.3.2 _ldType HeadStageIDType ^ _IdType
```

11.131.4 Property Documentation

```
\textbf{11.131.4.1} \quad \textbf{AdditionalText} \quad \texttt{String}^{\wedge} \; \texttt{AdditionalText} \quad \texttt{[get], [set]}
```

```
11.131.4.2 IdType HeadStageIDType^ IdType [get]
```

11.132 HeadStageIDTypeState Class Reference

Properties

- unsigned int State [get]
- HeadStageIDType [get]
- bool ControlState [get]
- bool DataState [get]

11.132.1 Property Documentation

```
11.132.1.1 ControlState bool ControlState [get]
```

11.132.1.2 DataState bool DataState [get]

11.132.1.3 IdType HeadStageIDType^ IdType [get]

11.132.1.4 State unsigned int State [get]

11.133 mkfilterNet Class Reference

Static Public Member Functions

- static int mkfilter (String[^] filtertype, double value, String[^] passtype, int order, double alpha1, double alpha2, [System::Runtime::InteropServices::Out] array< double >[^]% xcoeffs, [System::Runtime::InteropServices::Out] array< double >[^]% ycoeffs)
- static int mkfilter_MCS (int SamplesPerSecond, double R1, double R2, double C, double Amplification, double Correction, [System::Runtime::InteropServices::Out] array< double >^% xcoeffs, [System::Runtime::
 InteropServices::Out] array< double >^% ycoeffs)
- static int mkfilter_MCS (int SamplesPerSecond, double R1, double R2, double C, double Correction, [System::Runtime::InteropServices::Out] array< double >^% xcoeffs, [System::Runtime::InteropServices::Out] array< double >^% ycoeffs)
- static int mkfilter_MCS_k (int SamplesPerSecond, double R1, double R2, double C, double Amplification, double Correction, [System::Runtime::InteropServices::Out] array< double >^% coeffs)
- static int mkfilter_MCS_k (int SamplesPerSecond, double R1, double R2, double C, double Correction, [System::Runtime::InteropServices::Out] array< double >^% coeffs)
- static void mkfilter_coef_in_one_set (int n, [System::Runtime::InteropServices::In] array< double >^ xcoeffs, [System::Runtime::InteropServices::Out] array< double >^% out_coeffs)
- static void mkfilter_scale_coef_in_one_set (int n, double scale, [System::Runtime::InteropServices::In] array< double >^ xcoeffs, [System::Runtime::InteropServices::In] array< double >^ ycoeffs, [System::

 Runtime::InteropServices::Out] array< double >^% out_coeffs)
- static void mkfilter_normalize_coeffs_short (short maxvalue, [System::Runtime::InteropServices::In] array
 double >^ coeffs, [System::Runtime::InteropServices::Out] array< short >^% out_coeffs)
- static void mkfilter_normalize_coeffs_int (int maxvalue, [System::Runtime::InteropServices::In] array< double >^ coeffs, [System::Runtime::InteropServices::Out] array< int >^% out coeffs)
- static void mkfilter_normalize_scale_coeffs_int (int maxvalue, [System::Runtime::InteropServices::In] array
 double >^ coeffs, [System::Runtime::InteropServices::Out] array< int >^% out coeffs)
- static double mkfilter_highpass_coeff (int SamplesPerSecond, double Frequency)
- static double mkfilter_highpass_k (int SamplesPerSecond, double Frequency)
- static double mkfilter_highpass_frequency_from_coeff (int SamplesPerSecond, double coeff)
- static double mkfilter_highpass_frequency_from_k (int SamplesPerSecond, double k)

11.133.1 Member Function Documentation

```
11.133.1.2 mkfilter_coef_in_one_set() static void mkfilter_coef_in_one_set (
              int n_{i}
              [System::Runtime::InteropServices::In] array< double >^{\land} xcoeffs,
              [System::Runtime::InteropServices::In] array< double >^{\land} ycoeffs,
              [System::Runtime::InteropServices::Out] \ array<\ double > ^% \ out\_coeffs ) \ [static]
\textbf{11.133.1.3} \quad \textbf{mkfilter\_highpass\_coeff()} \quad \texttt{static double mkfilter\_highpass\_coeff} \ (
              int SamplesPerSecond,
              double Frequency ) [static]
11.133.1.4 mkfilter_highpass_frequency_from_coeff() static double mkfilter_highpass_frequency_←
from_coeff (
              int SamplesPerSecond,
              double coeff ) [static]
11.133.1.5 mkfilter_highpass_frequency_from_k() static double mkfilter_highpass_frequency_from←
_k (
              int SamplesPerSecond,
              double k ) [static]
11.133.1.6 mkfilter_highpass_k() static double mkfilter_highpass_k (
              int SamplesPerSecond,
              double Frequency ) [static]
11.133.1.7 mkfilter_MCS() [1/2] static int mkfilter_MCS (
              int SamplesPerSecond,
              double R1,
              double R2,
              double C_{\bullet}
              double Amplification,
              double Correction,
              [System::Runtime::InteropServices::Out] array< double >^% xcoeffs,
              [System::Runtime::InteropServices::Out] array<br/> double >^{\%} ycoeffs ) [static]
11.133.1.8 mkfilter_MCS() [2/2] static int mkfilter_MCS (
              int SamplesPerSecond,
              double R1,
              double R2,
              double C,
              double Correction,
              [System::Runtime::InteropServices::Out] array< double >^% xcoeffs,
              [System::Runtime::InteropServices::Out] array< double >^{\%} ycoeffs ) [static]
```

```
11.133.1.9 mkfilter_MCS_k() [1/2] static int mkfilter_MCS_k (
              int SamplesPerSecond,
             double R1,
             double R2,
             double C,
             double Amplification,
             double Correction,
             [System::Runtime::InteropServices::Out] array< double >^{\%} coeffs) [static]
11.133.1.10 mkfilter\_MCS\_k() [2/2] static int mkfilter\_MCS\_k (
             int SamplesPerSecond,
             double R1,
             double R2,
             double C_{\bullet}
             double Correction,
              [System::Runtime::InteropServices::Out] array< double >^{\%} coeffs ) [static]
11.133.1.11 mkfilter_normalize_coeffs_int() static void mkfilter_normalize_coeffs_int (
             int maxvalue,
             [System::Runtime::InteropServices::In] array< double >^{\wedge} coeffs,
              [System::Runtime::InteropServices::Out] array< int >^% out_coeffs ) [static]
11.133.1.12 mkfilter_normalize_coeffs_short() static void mkfilter_normalize_coeffs_short (
             short maxvalue.
              [System::Runtime::InteropServices::In] array< double >^{\wedge} coeffs,
              [System::Runtime::InteropServices::Out] array< short >^% out_coeffs) [static]
11.133.1.13 mkfilter_normalize scale coeffs_int() static void mkfilter_normalize_scale_coeffs_int
(
             int maxvalue,
              [System::Runtime::InteropServices::In] array< double >^{\land} coeffs,
              [System::Runtime::InteropServices::Out] array< int >^% out_coeffs ) [static]
11.133.1.14 mkfilter_scale_coef_in_one_set() static void mkfilter_scale_coef_in_one_set (
             int n_{i}
             double scale,
              [System::Runtime::InteropServices::In] array< double >^{\land} xcoeffs,
              [System::Runtime::InteropServices::In] array< double >^{\land} ycoeffs,
              [System::Runtime::InteropServices::Out] \ array<\ double > ^\% \ out\_coeffs \ ) \ [static]
```

11.134 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference

Public Member Functions

- · void SetParameter (unsigned short command, unsigned short index, unsigned int value)
- · void SetParameter (unsigned short command, unsigned short index, unsigned int value1, unsigned int value2)
- void SetUserParameter (unsigned short index, unsigned int value)

Stores persistently 32 bit integer values on RoboMain

void SetUserParameter (unsigned short index, int value)

Stores persistently 32 bit integer values on RoboMain

- void GetParameter (unsigned short command, unsigned short index, [System::Runtime::InteropServices::
 Out]unsigned int% value)
- void GetParameter (unsigned short command, unsigned short index, [System::Runtime::InteropServices::
 Out]unsigned int% value1, [System::Runtime::InteropServices::Out]unsigned int% value2)
- $\bullet \ \ void\ Get User Parameter\ (unsigned\ short\ index,\ [System::Runtime::InteropServices::Out] unsigned\ int\%\ value)$

Reads 32 bit integer values stored persistently on RoboMain

void GetUserParameter (unsigned short index, [System::Runtime::InteropServices::Out]int% value)

Reads 32 bit integer values stored persistently on RoboMain

- void FindReferencePhase0 (unsigned char busaddress, char axes)
- void FindReferencePhase0 (unsigned char busaddress, char axes, int timeout)
- unsigned char HasRef (unsigned char busaddress, char axes)
- void SetHWRevision (unsigned int revision)
- unsigned int GetHWRevision ()
- void SetHWConfig (unsigned int config)
- unsigned int GetHWConfig ()
- void SetMinPressureWaitTime (unsigned int t)
- unsigned int GetMinPressureWaitTime ()
- void SetMinPressure (unsigned int pressure)
- unsigned int GetMinPressure ()
- void SetMaxPressureWaitTime (unsigned int t)
- unsigned int GetMaxPressureWaitTime ()
- void SetMinNoPressureWaitTime (unsigned int t)
- unsigned int GetMinNoPressureWaitTime ()
- void SetMaxNoPressure (unsigned int pressure)
- unsigned int GetMaxNoPressure ()
- void SetMaxNoPressureWaitTime (unsigned int t)
- unsigned int GetMaxNoPressureWaitTime ()
- · void SetSearchReferenceMethod (unsigned char busaddress, char axes, unsigned int method)
- unsigned int GetSearchReferenceMethod (unsigned char busaddress, char axes)
- void SetSearchReferenceOffsetPos (unsigned char busaddress, char axes, int offsetpos)
- int GetSearchReferenceOffsetPos (unsigned char busaddress, char axes)
- void SetSearchReferenceFastSpeed (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short GetSearchReferenceFastSpeed (unsigned char busaddress, char axes)
- void SetSearchReferenceFastAccel (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short GetSearchReferenceFastAccel (unsigned char busaddress, char axes)
- void SetSearchReferenceFineSpeed (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short GetSearchReferenceFineSpeed (unsigned char busaddress, char axes)
- void SetSearchReferenceFineAccel (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short GetSearchReferenceFineAccel (unsigned char busaddress, char axes)
- void SetSearchReferenceMoveOut (unsigned char busaddress, char axes, int move)
- int GetSearchReferenceMoveOut (unsigned char busaddress, char axes)
- void SetAxisConfig (unsigned char busaddress, char axes, unsigned int config)
- unsigned int GetAxisConfig (unsigned char busaddress, char axes)
- void GetPhases (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out] unsigned short% phase0, [System::Runtime::InteropServices::Out] unsigned short% lastphase)

11.134.1 Member Function Documentation

```
11.134.1.1 FindReferencePhaseO() [1/2] void FindReferencePhaseO (
            unsigned char busaddress,
            char axes )
11.134.1.2 FindReferencePhaseO() [2/2] void FindReferencePhaseO (
            unsigned char busaddress,
            char axes,
            int timeout )
11.134.1.3 GetAxisConfig() unsigned int GetAxisConfig (
            unsigned char busaddress,
             char axes )
11.134.1.4 GetHWConfig() unsigned int GetHWConfig ( )
11.134.1.5 GetHWRevision() unsigned int GetHWRevision ()
11.134.1.6 GetMaxNoPressure() unsigned int GetMaxNoPressure ( )
11.134.1.7 GetMaxNoPressureWaitTime() unsigned int GetMaxNoPressureWaitTime ( )
11.134.1.8 GetMaxPressureWaitTime() unsigned int GetMaxPressureWaitTime ( )
11.134.1.9 GetMinNoPressureWaitTime() unsigned int GetMinNoPressureWaitTime ()
```

```
11.134.1.10 GetMinPressure() unsigned int GetMinPressure ( )
11.134.1.11 GetMinPressureWaitTime() unsigned int GetMinPressureWaitTime ( )
11.134.1.12 GetParameter() [1/2] void GetParameter (
             unsigned short command,
             unsigned short index,
             [{\tt System::Runtime::InteropServices::Out}] \ {\tt unsigned int \$} \ {\tt value} \ )
11.134.1.13 GetParameter() [2/2] void GetParameter (
             unsigned short command,
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned int% value1,
             [System::Runtime::InteropServices::Out] unsigned int% value2 )
11.134.1.14 GetPhases() void GetPhases (
             unsigned char busaddress,
             char axes,
             [System::Runtime::InteropServices::Out] unsigned short% phase0,
             [System::Runtime::InteropServices::Out] unsigned short% lastphase )
11.134.1.15 GetSearchReferenceFastAccel() unsigned short GetSearchReferenceFastAccel (
             unsigned char busaddress,
             char axes )
11.134.1.16 GetSearchReferenceFastSpeed() unsigned short GetSearchReferenceFastSpeed (
             unsigned char busaddress,
             char axes )
11.134.1.17 GetSearchReferenceFineAccel() unsigned short GetSearchReferenceFineAccel (
             unsigned char busaddress,
             char axes )
```

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

index	address offset of parameter; range: 015
value	data buffer

```
11.134.1.23 GetUserParameter() [2/2] void GetUserParameter (
unsigned short index,

[System::Runtime::InteropServices::Out] unsigned int% value )
```

[System::Runtime::InteropServices::Out] int% value)

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

index	address offset of parameter; range: 015
value	data buffer

```
11.134.1.24 HasRef() unsigned char HasRef (
            unsigned char busaddress,
             char axes )
11.134.1.25 SetAxisConfig() void SetAxisConfig (
            unsigned char busaddress,
            char axes,
             unsigned int config )
11.134.1.26 SetHWConfig() void SetHWConfig (
             unsigned int config )
11.134.1.27 SetHWRevision() void SetHWRevision (
            unsigned int revision )
11.134.1.28 SetMaxNoPressure() void SetMaxNoPressure (
            unsigned int pressure )
11.134.1.29 SetMaxNoPressureWaitTime() void SetMaxNoPressureWaitTime (
            unsigned int t )
11.134.1.30 SetMaxPressureWaitTime() void SetMaxPressureWaitTime (
             unsigned int t )
11.134.1.31 SetMinNoPressureWaitTime() void SetMinNoPressureWaitTime (
            unsigned int t )
11.134.1.32 SetMinPressure() void SetMinPressure (
             unsigned int pressure )
```

```
11.134.1.33 SetMinPressureWaitTime() void SetMinPressureWaitTime (
             unsigned int t )
11.134.1.34 SetParameter() [1/2] void SetParameter (
             unsigned short command,
             unsigned short index,
             unsigned int value )
11.134.1.35 SetParameter() [2/2] void SetParameter (
            unsigned short command,
             unsigned short index,
             unsigned int value1,
             unsigned int value2)
11.134.1.36 SetSearchReferenceFastAccel() void SetSearchReferenceFastAccel (
             unsigned char busaddress,
             char axes,
             unsigned short accel )
11.134.1.37 SetSearchReferenceFastSpeed() void SetSearchReferenceFastSpeed (
             unsigned char busaddress,
             char axes,
             unsigned short speed )
11.134.1.38 SetSearchReferenceFineAccel() void SetSearchReferenceFineAccel (
             unsigned char busaddress,
             char axes,
             unsigned short accel )
11.134.1.39 SetSearchReferenceFineSpeed() void SetSearchReferenceFineSpeed (
             unsigned char busaddress,
             char axes,
             unsigned short speed )
```

$\textbf{11.134.1.40} \quad \textbf{SetSearchReferenceMethod()} \quad \texttt{void SetSearchReferenceMethod} \quad \textbf{(}$

```
unsigned char busaddress,
char axes,
unsigned int method )
```

$\textbf{11.134.1.41} \quad \textbf{SetSearchReferenceMoveOut()} \quad \texttt{void SetSearchReferenceMoveOut} \quad \textbf{(}$

```
unsigned char busaddress,
char axes,
int move )
```

11.134.1.42 SetSearchReferenceOffsetPos() void SetSearchReferenceOffsetPos (

```
unsigned char busaddress,
char axes,
int offsetpos )
```

11.134.1.43 SetUserParameter() [1/2] void SetUserParameter (

```
unsigned short index,
int value )
```

Stores persistently 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

index address offset of parame		address offset of parameter; range: 015
	value	data to be stored

$\textbf{11.134.1.44} \quad \textbf{SetUserParameter() [2/2]} \quad \texttt{void SetUserParameter ()}$

```
unsigned short index,
unsigned int value )
```

Stores persistently 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

index	(address offset of parameter; range: 015
value	,	data to be stored

11.135 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference

Public Member Functions

- void FindReferencePhase0XY ()
- void FindReferencePhase0XY (int timeout)

11.135.1 Member Function Documentation

```
11.135.1.1 FindReferencePhaseOXY() [1/2] void FindReferencePhaseOXY ( )
```

```
11.135.1.2 FindReferencePhaseOXY() [2/2] void FindReferencePhaseOXY ( int timeout )
```

11.136 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference

Public Member Functions

- s_FilterAttributesNet (s_FilterAttributes attrib)
- s_FilterAttributes ToCpp ()

Public Attributes

- uint32_t PreCommaB
- uint32_t PostCommaB
- uint32_t CommaPositionB
- uint32_t PreCommaA
- uint32_t PostCommaA
- uint32_t CommaPositionA

11.136.1 Constructor & Destructor Documentation

11.136.2 Member Function Documentation

11.136.2.1 ToCpp() s_FilterAttributes ToCpp ()

11.136.3 Member Data Documentation

11.136.3.1 CommaPositionA uint32_t CommaPositionA

11.136.3.2 CommaPositionB uint32_t CommaPositionB

11.136.3.3 PostCommaA uint32_t PostCommaA

11.136.3.4 PostCommaB uint32_t PostCommaB

11.136.3.5 PreCommaA uint32_t PreCommaA

11.136.3.6 PreCommaB uint32_t PreCommaB

11.137 CMeaAudioFunctionNet::s_setaudionet Struct Reference

Public Attributes

- int channel
- · int amplification

11.137.1 Member Data Documentation

$\textbf{11.137.1.1} \quad \textbf{amplification} \quad \text{int amplification} \\$

```
11.137.1.2 channel int channel
```

11.138 CStimulusFunctionNet::SidebandData Class Reference

Public Member Functions

- SidebandData ()
- ∼SidebandData ()

Destructor: called by Dispose()

• !SidebandData ()

Finalizer: called by GC before collecting

Properties

```
• array< int32_t >^{\land} Sideband [get]
```

• array< uint64_t >^ Duration [get]

11.138.1 Constructor & Destructor Documentation

```
11.138.1.1 SidebandData() SidebandData ()
```

```
11.138.1.2 \simSidebandData() \simSidebandData ()
```

Destructor: called by Dispose()

```
11.138.1.3 "!SidebandData() !SidebandData ( )
```

Finalizer: called by GC before collecting

11.138.2 Property Documentation

```
11.138.2.1 Duration array< uint64_t>^{\wedge} Duration [get]
```

11.138.2.2 Sideband array< int32_t> $^{\land}$ Sideband [get]

11.139 StgStatusNet Class Reference

Static Public Member Functions

- static StgStatusNet ^ FromIntPtr (IntPtr stgstatus)
- static StgStatusNet ^ FromPtr (stgstatus_t *stgstatus)

Public Attributes

- array< Stg200xTriggerStatusEnumNet > ^ TiggerStatus
- array< uint32_t > $^{\land}$ ListOfChangedTriggers

11.139.1 Member Function Documentation

11.139.2 Member Data Documentation

11.139.2.1 ListOfChangedTriggers array<uint32_t> ^ ListOfChangedTriggers

 $\textbf{11.139.2.2} \quad \textbf{TiggerStatus} \quad \texttt{array} < \texttt{Stg200xTriggerStatusEnumNet} > \\ \land \quad \texttt{TiggerStatus}$

11.140 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference

Public Member Functions

- StimulusDeviceDataAndUnrolledData ()
- $\bullet \ \sim \! Stimulus Device Data And Unrolled Data ()$

Destructor: called by Dispose()

• !StimulusDeviceDataAndUnrolledData ()

Finalizer: called by GC before collecting

Properties

```
• array< uint8_t >^ DeviceData [get]
   • int DeviceDataLength [get]
    • array< int32_t >^ UnrolledAmplitude [get]
    • array< uint32 t >^{\wedge} UnrolledSync [get]
    • array< uint64_t >^ UnrolledDuration [get]
11.140.1 Constructor & Destructor Documentation
11.140.1.1 StimulusDeviceDataAndUnrolledData() StimulusDeviceDataAndUnrolledData ()
11.140.1.2 ~StimulusDeviceDataAndUnrolledData() ~StimulusDeviceDataAndUnrolledData ()
Destructor: called by Dispose()
11.140.1.3 "!StimulusDeviceDataAndUnrolledData() !StimulusDeviceDataAndUnrolledData ()
Finalizer: called by GC before collecting
11.140.2 Property Documentation
11.140.2.1 DeviceData array< uint8_{t}>^ DeviceData [get]
11.140.2.2 DeviceDataLength int DeviceDataLength [get]
11.140.2.3 UnrolledAmplitude array< int32_t>^ UnrolledAmplitude [get]
\textbf{11.140.2.4} \quad \textbf{UnrolledDuration} \quad \texttt{array} < \ \texttt{uint64\_t} > ^{\wedge} \ \texttt{UnrolledDuration} \quad \texttt{[get]}
```

11.140.2.5 UnrolledSync array< $uint32_t>^{\wedge} UnrolledSync$ [get]

11.141 usbSetupPacket_t Class Reference

Public Attributes

- uint8_t bmRequestType
- uint8_t bRequest
- uint16_t wValue
- uint16_t wlndex
- · uint16_t wLength

11.141.1 Member Data Documentation

```
11.141.1.1 bmRequestType uint8_t bmRequestType
```

```
11.141.1.2 bRequest uint8_t bRequest
```

```
11.141.1.3 windex uint16_t wIndex
```

```
11.141.1.4 wLength uint16_t wLength
```

11.141.1.5 wValue uint16_t wValue

11.142 W2100_StimulusParametersNet Struct Reference

Public Attributes

- int DACResolution
- int TimeResolutionInNanoSeconds
- int VoltageRangeInMicroVolt
- int VoltageResolutionInMicroVolt
- int CurrentRangeInNanoAmp
- int CurrentResolutionInNanoAmp

11 142 1	Member	Data Do	cumen	tation

- 11.142.1.1 CurrentRangeInNanoAmp int CurrentRangeInNanoAmp
- 11.142.1.2 CurrentResolutionInNanoAmp int CurrentResolutionInNanoAmp
- 11.142.1.3 DACResolution int DACResolution
- 11.142.1.4 TimeResolutionInNanoSeconds int TimeResolutionInNanoSeconds
- $\textbf{11.142.1.5} \quad \textbf{VoltageRangeInMicroVolt} \quad \texttt{int VoltageRangeInMicroVolt}$
- $\textbf{11.142.1.6} \quad \textbf{VoltageResolutionInMicroVolt} \quad \texttt{int VoltageResolutionInMicroVolt}$

Index

!CDacCalibrationFunctionNet	CWarnerUssingDeviceNet, 540
CDacCalibrationFunctionNet, 51	!CWarnerUssingFunctionNet
!CDigOutStimulatorFunctionNet	CWarnerUssingFunctionNet, 543
CDigOutStimulatorFunctionNet, 60	!CWarnerValveControllerDeviceNet
!CExternDTesterDeviceNet	CWarnerValveControllerDeviceNet, 560
CExternDTesterDeviceNet, 65	!CWarnerValveControllerDeviceTesterFunctionNet
!CGrapheneFunctionNet	CWarnerValveControllerDeviceTesterFunctionNet,
CGrapheneFunctionNet, 103	581
!CInterfaceboard2FunctionNet	!SidebandData
CInterfaceboard2FunctionNet, 117	CStimulusFunctionNet::SidebandData, 617
!CInterfaceboardFunctionNet	!StimulusDeviceDataAndUnrolledData
CInterfaceboardFunctionNet, 119	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
!CLIH3DeviceNet	619
CLIH3DeviceNet, 122	AdditionalText
!CMEA2100x256FunctionNet	HeadstageIDTypeObject, 603
CMEA2100x256FunctionNet, 262	_ldType
!CMcsUsbFunctionNet	HeadstageIDTypeObject, 604
CMcsUsbFunctionNet, 231	~CCMOSMeaDeviceNet
!CMcsUsbListNet	CCMOSMeaDeviceNet, 46
CMcsUsbListNet, 238	~CChannelTestDeviceNet
!CMcsUsbNet	CChannelTestDeviceNet, 33
CMcsUsbNet, 244	~CCreateFilterNet
!CMeFunctionNet	CCreateFilterNet, 48
CMeFunctionNet, 296	~CDacCalibrationFunctionNet
!CMeaCleanDeviceNet	CDacCalibrationFunctionNet, 51
CMeaCleanDeviceNet, 268	~CDigOutStimulatorFunctionNet
!CMeaCoatDeviceNet	CDigOutStimulatorFunctionNet, 60
	~CExternDTesterDeviceNet
CMeaCoatDeviceNet, 272	
!CMultiBatteryChargerDeviceNet	CExternDTesterDeviceNet, 65
CMultiBatteryChargerDeviceNet, 298	~CFilterCoefficientsNet
!CMultiwellCallbackFunctionNet	CFilterCoefficientsNet, 67
CMultiwellCallbackFunctionNet, 305	~CFilterPropertyNet
!CMultiwellDeviceNet	CFilterPropertyNet, 72
CMultiwellDeviceNet, 307	~CFluidControlDeviceNet
!CMultiwellOptoStimFunctionNet	CFluidControlDeviceNet, 75
CMultiwellOptoStimFunctionNet, 313	~CGenericDevelopDeviceNet
!CPPCFunctionNet	CGenericDevelopDeviceNet, 89
CPPCFunctionNet, 347	~CGilsonDeviceNet
!CPedoterDeviceNet	CGilsonDeviceNet, 101
CPedoterDeviceNet, 328	~CGrapheneFunctionNet
!CPositionIIDeviceNet	CGrapheneFunctionNet, 103
CPositionIIDeviceNet, 334	\sim CInterfaceboard2FunctionNet
!CPositionImpDeviceNet	CInterfaceboard2FunctionNet, 117
CPositionImpDeviceNet, 342	\sim CInterfaceboardFunctionNet
!CProgramPressureCurveNet	CInterfaceboardFunctionNet, 119
CProgramPressureCurveNet, 361	\sim CLIH3DeviceNet
!CPulseGeneratorFunctionNet	CLIH3DeviceNet, 122
CPulseGeneratorFunctionNet, 362	\sim CMEA2100x256FunctionNet
!CRFFunctionNet	CMEA2100x256FunctionNet, 262
CRFFunctionNet, 370	\sim CMcsBusNet
!CSCUFunctionNet	CMcsBusNet, 168
CSCUFunctionNet, 423	\sim CMcsBus_AxisParametersNet
!CTEERFunctionNet	CMcsBus_AxisParametersNet, 130
CTEERFunctionNet, 513	~CMcsBus_ExtensionNet
!CWarnerUssingDeviceNet	CMcsBus ExtensionNet. 131

~CMcsBus_FYIExtensionNet	~CPositionIIDeviceNet
CMcsBus_FYIExtensionNet, 132	CPositionIIDeviceNet, 334
~CMcsBus_MotorControlNet	\sim CPositionImpDeviceNet
CMcsBus_MotorControlNet, 136	CPositionImpDeviceNet, 342
~CMcsBus_SensorNet	\sim CProgramPressureCurveNet
CMcsBus_SensorNet, 152	CProgramPressureCurveNet, 360
~CMcsBus_TempSensorNet	~CPulseGeneratorFunctionNet
CMcsBus_TempSensorNet, 161	CPulseGeneratorFunctionNet, 362
~CMcsBus VoltageModeNet	~CRFFunctionNet
CMcsBus VoltageModeNet, 164	CRFFunctionNet, 370
~CMcsUsbDacqNet	~CRetinaLedDeviceNet
CMcsUsbDacqNet, 177	CRetinaLedDeviceNet, 368
~CMcsUsbFactoryNet	~CRoboDeviceNet
CMcsUsbFactoryNet, 223	CRoboDeviceNet, 394
~CMcsUsbFunctionNet	~CRoboFluidDeviceNet
CMcsUsbFunctionNet, 231	CRoboFluidDeviceNet, 406
~CMcsUsbListEntryNet	~CSCUFunctionNet
•	
CMcsUsbListEntryNet, 232 ~CMcsUsbListNet	CSCUFunctionNet, 423
	~CSafeISDeviceNet
CMcsUsbListNet, 238	CSafeISDeviceNet, 418
~CMcsUsbNet	~CStg200xBasicNet
CMcsUsbNet, 244	CStg200xBasicNet, 442
~CMeFunctionNet	~CStg200xDownloadNet
CMeFunctionNet, 296	CStg200xDownloadNet, 478
~CMeaCleanDeviceNet	~CSw2to64DeviceNet
CMeaCleanDeviceNet, 267	CSw2to64DeviceNet, 496
~CMeaCoatDeviceNet	\sim CTEERFunctionNet
CMeaCoatDeviceNet, 272	CTEERFunctionNet, 513
\sim CMeaDeviceNet	~CTEERMachineDeviceNet
CMeaDeviceNet, 278	CTEERMachineDeviceNet, 522
\sim CMeaImpedanceDeviceNet	\sim CTcxDeviceNet
CMealmpedanceDeviceNet, 290	CTcxDeviceNet, 500
\sim CMeaSwitchDeviceNet	~CWarnerUssingDeviceNet
CMeaSwitchDeviceNet, 293	CWarnerUssingDeviceNet, 540
\sim CMeaUSBDeviceNet	~CWarnerUssingFunctionNet
CMeaUSBDeviceNet, 295	CWarnerUssingFunctionNet, 543
~CMultiBatteryChargerDeviceNet	~CWarnerValveControllerDeviceNet
CMultiBatteryChargerDeviceNet, 298	CWarnerValveControllerDeviceNet, 560
~CMultiwellCallbackFunctionNet	~CWarnerValveControllerDeviceTesterFunctionNet
CMultiwellCallbackFunctionNet, 305	CWarnerValveControllerDeviceTesterFunctionNet,
~CMultiwellDeviceNet	581
CMultiwellDeviceNet, 307	~DriverVersionNet
~CMultiwellOptoStimFunctionNet	DriverVersionNet, 592
CMultiwellOptoStimFunctionNet, 313	~SidebandData
~CNF_GenDeviceNet	CStimulusFunctionNet::SidebandData, 617
CNF_GenDeviceNet, 317	~StimulusDeviceDataAndUnrolledData
~COkuvisionStimulatorDeviceNet	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData.
COkuvisionStimulatorDeviceNet, 323	619
~CPPCFunctionNet	010
CPPCFunctionNet, 347	A
~CPathIdentDeviceNet	CFilterCoefficientsNet, 68
	AdditionalText
CPathIdentDeviceNet, 327 ~CPedoterDeviceNet	HeadstageIDTypeObject, 604
	AddLoopEntry
CPerioteltic Pump Povice Net	CRetinaLedDeviceNet, 368
~CPeristalticPumpDeviceNet	AddSelectedChannelsQueue
CPeristalticPumpDeviceNet, 329	CMcsUsbDacqNet, 177, 179, 180
~CPgaDeviceNet	AddSoftwareKey
CPgaDeviceNet, 330	CMcslishNat 244

AddTableEntry	usbSetupPacket_t, 620
CRetinaLedDeviceNet, 368	BurnAdcOffset
Altera	COctoPotDeviceNet, 319
FirmwareDestinationNames, 597	BurnDacOffset
amplification	CDacCalibrationFunctionNet, 51
CMeaAudioFunctionNet::s_setaudionet, 616	COctoPotDeviceNet, 319
CW2100_FunctionNet::AudioChannelsNet, 29	BUS1_MCSBUS1
AmplifierSettle	FirmwareDestinationNames, 597
CIntanMea_FunctionNet, 115	BUS1_MCSBUS2
AnalogGain	FirmwareDestinationNames, 597
CMeaDeviceNet, 283	BusType
ApplyGains	DeviceIdNet, 590
CPgaDeviceNet, 331	ButterworthFilterHighPassNet, 31
AreTransistorVoltagesSet	ButterworthFilterHighPassNet, 32
CCMOSMea_FunctionNet, 36	ButterworthFilterLowPassNet, 32
AssociateToThis	ButterworthFilterLowPassNet, 32
CMcsUsbNet, 244	0.111
AutomaticAnalogOut	CalibrateThermocouple
CSCUFunctionNet, 423	CFluidControlDeviceNet, 75
Axes_I	CTcxDeviceNet, 500
CRoboDeviceNet, 401	CancelInternalCalibration
Axes_X	CTEERFunctionNet, 513
CRoboDeviceNet, 401	CancelPoolLoop
Axes_Y	CRoboDeviceNet, 394
CRoboDeviceNet, 401	CancelPoolLoopAndStopMovement
Axes_Z	CRoboDeviceNet, 394
CRoboDeviceNet, 401	CancelTableLoop
Axis_I	CRoboDacqNet, 380
CRoboDeviceNet, 401	CancelTableLoopAndStopTable
Axis_X	CRoboDacqNet, 380
CRoboDeviceNet, 401	CapacityTest
Axis_Y	CMultiBatteryChargerDeviceNet, 298 CatchAmpGetAdcMean
CRoboDeviceNet, 401	CMcsBus SensorNet, 152
Axis_Z	CatchAmpGetAdcValue
CRoboDeviceNet, 401	CMcsBus_SensorNet, 153
В	CatchAmpGetAdcValueH
CFilterCoefficientsNet, 68	CMcsBus_SensorNet, 153
BatteryState, 29	CatchAmpGetAdcValueL
Charge, 30	CMcsBus SensorNet, 153
ChargeRegionString, 30	CatchAmpGetDacAmplitude
ChargeString, 30	CMcsBus_SensorNet, 153
Voltage, 30	CatchAmpGetDacEnable
VoltageString, 30	CMcsBus_SensorNet, 153
BcdDevice	CatchAmpGetDacOffset
DeviceIdNet, 590	CMcsBus_SensorNet, 153
BeginImpedanceCheck	CatchAmpGetPwmEnable
CIntanMea_FunctionNet, 115	CMcsBus_SensorNet, 153
BesselFilterHighPassNet, 30	CatchAmpSetDacAmplitude
BesselFilterHighPassNet, 30	CMcsBus_SensorNet, 153
BesselFilterLowPassNet, 31	CatchAmpSetDacEnable
BesselFilterLowPassNet, 31	CMcsBus_SensorNet, 154
bmRequestType	CatchAmpSetDacOffset
usbSetupPacket_t, 620	CMcsBus_SensorNet, 154
BOOST_BIT	CatchAmpSetPwmEnable
CW2100_StimulatorFunctionNet, 538	CMcsBus_SensorNet, 154
Bootstrap	CChannelTestDeviceNet, 33
FirmwareDestinationNames, 597	~CChannelTestDeviceNet, 33
bRequest	CChannelTestDeviceNet, 33

Cot Amplitudo 22	UndataTransistarValtages 44
SetAmplitude, 33	UpdateTransistorVoltages, 44
SetAttenuation, 33	VOPSTimerSetResetTimes, 44, 45
SetFrequency, 33	CCMOSMeaDeviceNet, 45
SetWaveform, 34	~CCMOSMeaDeviceNet, 46
CCMOSMea_FunctionNet, 34	CCMOSMeaDeviceNet, 46
AreTransistorVoltagesSet, 36	CMosMea, 47
CCMOSMea_FunctionNet, 36	GetAvailableBaseSamplerates, 46
ClearSTGOutput, 36	GetBaseSamplerate, 46
DetectChipType, 36	GetChannelDatal16, 46
EnableChannelsInGroup, 36, 37	GetChannelDatal32, 46
GetADCInputOffset, 37	GetChannelDataUI16, 46
GetBath, 37	GetChannelDataUI32, 47
GetBathMode, 37	GetCMOSDataDictionary, 47
GetEnabledChannelsInGroup, 37	SetBaseSamplerate, 47
GetGate, 37	SetRegionOfInterests, 47
GetGNDI, 37	Stimulus, 47
GetGroupADCBits, 37	UpdateChannelBlock, 47
GetGroupChannelBitmaskBySelect, 38	CCMOSMeaDeviceNet::CRegionOfInterestRect, 366
GetGroupChannelBitmaskHS1NCBathCurrent, 38	CRegionOfInterestRect, 366
GetGroupChannelBitmaskHS1NCCol2Current, 38	DeepCopy, 366
GetGroupChannelBitmaskHS1NChipTemp, 38	m_Bottom, 367
GetGroupChannelBitmaskHS1Sidebands, 39	m_Left, 367
GetGroupChannelBitmaskHS1TriggerStatus, 39	m_Right, 367
GetGroupChannelBitmaskIFDigChannels, 39	m_Top, 367
GetGroupChannelBitmaskInterfaceADC, 39	CCreateFilterNet, 48
GetGroupChannelBitmaskPacketFrameContext, 40	\sim CCreateFilterNet, 48
GetGroupChannelBitmaskSTG1DACSignal, 40	CCreateFilterNet, 48, 49
GetGroupDCOffset, 40	CutoffFrequency, 49
GetGroupID, 40	FindFilter, 49
GetGroupNumberOfChannels, 41	GetBiQuad, 49
GetGroupResolutionPerDigit, 41	GetBiQuads, 49
GetGroupSampleSize, 41	NumCoefSets, 49
GetGroupType, 41	Order, 49
GetGroupUnit, 41, 42	SampleRate, 50
GetMaxNumOfColumns, 42	Scale, 50
GetNeurochipMemoryData, 42	CDacCalibrationFunctionNet, 50
GetNeurochipMemorySize, 42	!CDacCalibrationFunctionNet, 51
GetNumberOfSupportedGroups, 42	~CDacCalibrationFunctionNet, 51
GetSourceBulk, 42	BurnDacOffset, 51
GetSourceDrain, 42	CDacCalibrationFunctionNet, 50, 51
GetSourceGate, 42	GetDacOffset, 51
GetStimulusSites, 43	SetDacOffset, 51
GetVDD3I, 43	CDacqGroupChannelGenericSelectionNet, 52
GetVDDI, 43	CDacqGroupChannelGenericSelectionNet, 52
IsChipPowered, 43	CDacqGroupChannelSelectionNet, 52
IsGateFloating, 43	CDacqGroupChannelSelectionNet, 53
PowerChip, 43	CDacqGroupChannelSelectionTemplateNet
SetADCInputOffset, 43	CDacqGroupChannelSelectionTemplateNet< Dac-
•	
SetBath, 43	qGroupChannelEnumTemplateNet, Dac-
SetBathMode, 43	qGroupChannelEnumTemplate, CDevice-
SetGate, 43	GroupChannelInfoTemplateNet >, 54
SetGateFloating, 44	CDacqGroupChannelSelectionTemplateNet< Dac-
SetGateToVOP, 44	qGroupChannelEnumTemplateNet, Dac-
SetNeurochipMemoryData, 44	qGroupChannelEnumTemplate, CDevice-
SetSourceBulk, 44	GroupChannelInfoTemplateNet >, 53
SetSourceDrain, 44	CDacqGroupChannelSelectionTemplateNet, 54
SetSourceGate, 44	EnableChannelsInGroup, 54
SetStimulusSites, 44	GetDeviceGroupChannelInfos, 54

GetEnabledChannelsInGroup, 54	CommaPositionB, 616
GetGroupID, 54, 55	PostCommaA, 616
GetGroupNumberOfChannels, 55	PostCommaB, 616
GetGroupSampleSize, 55	PreCommaA, 616
GetGroupType, 55	PreCommaB, 616
GetNumberOfSupportedGroups, 55	s_FilterAttributesNet, 615
CDeviceGroupChannelInfoGenericNet, 56	ToCpp, 615
CDeviceGroupChannelInfoGenericNet, 56	CFilterConfigurationNet, 68
CDeviceGroupChannelInfoNet, 56	CFilterConfigurationNet, 68
CDeviceGroupChannelInfoNet, 57	EraseFilterParameterPermanent, 69
CDeviceGroupChannelInfoSCUNet, 57	GetFilterAttributes, 69
CDeviceGroupChannelInfoSCUNet, 57	GetHighpassFilterEnable, 69
CDeviceGroupChannelInfoTemplateNet	ResetHighpassFilter, 69
CDeviceGroupChannelInfoTemplateNet< Dacq-	SetFilterParameter, 69
GroupChannelEnumTemplateNet $>$, 58	SetFilterParameterPermanent, 69
CDeviceGroupChannelInfoTemplateNet< Dacq-	SetHighpassFilterEnable, 70
GroupChannelEnumTemplateNet $>$, 57	CFilterConfigurationRegisterNet, 70
CDeviceGroupChannelInfoTemplateNet, 58	CFilterConfigurationRegisterNet, 70
GroupID, 58	EraseFilterParameterPermanent, 70, 71
GroupType, 58	SetFilterParameter, 71
NumberOfChannels, 58	SetFilterParameterPermanent, 71
CDeviceGroupChannelInfoW2100Net, 58	CFilterPropertyNet, 72
CDeviceGroupChannelInfoW2100Net, 59	∼CFilterPropertyNet, 72
CDigOutStimulatorFunctionNet, 59	CFilterPropertyNet, 72
!CDigOutStimulatorFunctionNet, 60	CornerFrequency, 72
~CDigOutStimulatorFunctionNet, 60	CornerFrequencymHz, 73
CDigOutStimulatorFunctionNet, 60	FilterActive, 73
ClearChannel, 60	FilterBand, 73
GetGlobalRepeat, 60	FilterFamily, 73
GetNumberOfChannels, 61	FilterType, 73
GetStartTriggerSlope, 61	Order, 73
GetStopTriggerSlope, 61	ToString, 72
PrepareChannelData, 62	CFluidControlDeviceNet, 73
SendChannelData, 62	~CFluidControlDeviceNet, 75
SetGlobalRepeat, 62	CalibrateThermocouple, 75
SetStartTriggerSlope, 62	CFluidControlDeviceNet, 75
SetStopTriggerSlope, 63	GetAdc, 75
CEncapsulatorDeviceNet, 63	GetDigin, 75
CEncapsulatorDeviceNet, 64	GetDigout, 76
GetRoboFluidDevice, 64	GetPWM, 76
CExternDTesterDeviceNet, 64	GetReferenceTemperature, 76
!CExternDTesterDeviceNet, 65	GetSingleValve, 76
~CExternDTesterDeviceNet, 65	GetThermocoupleCalibration, 77
CExternDTesterDeviceNet, 65	GetThermocoupleNanovoltPerKelvin, 77
Read, 65	GetThermocoupleTemperature, 77
Read2, 65	GetValve, 78
Write, 65	McsBus_VoltageMode, 80
Write2, 66	SetDigout, 78
CFilterCoefficientsNet, 66	SetPWM, 78
∼CFilterCoefficientsNet, 67	SetSingleValve, 78
A, 68	SetThermocoupleNanovoltPerKelvin, 80
B, 68	SetValve, 80
CFilterCoefficientsNet, 66, 67	CFYIDeviceNet, 80
GetUintA, 67	CFYIDeviceNet, 81
GetUintB, 67	FYIProgram, 81
IsEqual, 67	FYITemp, 81
CFilterCoefficientsNet::s_FilterAttributesNet, 615	Sensor, 81
CommaPositionA, 616	CGenericDevelopDeviceNet, 82

- CGanaria Davialan Daviac Not. 90	CMaallahDaagNat 192 193
~CGenericDevelopDeviceNet, 89	CMcsUsbDacqNet, 182, 183
CGenericDevelopDeviceNet, 89	ChannelBlock_ReadAsFrameArrayUI16
ClosePipe, 89	CMcsUsbDacqNet, 183, 184
GetBuffer, 89	ChannelBlock_ReadAsFrameArrayUI32
GetByteBuffer, 89	CMcsUsbDacqNet, 185, 186
GetIntBuffer, 90	ChannelBlock_ReadFramesDictI16
GetShortBuffer, 91	CMcsUsbDacqNet, 186
GetUByteBuffer, 91	ChannelBlock ReadFramesDictl32
GetUIntBuffer, 93	CMcsUsbDacqNet, 187
GetUShortBuffer, 94	ChannelBlock ReadFramesDictUI16
OpenPipe, 94	CMcsUsbDacqNet, 187
ReadPipe, 95	ChannelBlock ReadFramesDictUI32
ResetPipe, 95	CMcsUsbDacqNet, 188
SetBuffer, 96	ChannelBlock ReadFramesI16
	-
SetByteBuffer, 96	Change Plant, Pand France 199
SetIntBuffer, 96	ChannelBlock_ReadFramesI32
SetShortBuffer, 97	CMcsUsbDacqNet, 190, 191
SetUByteBuffer, 97	ChannelBlock_ReadFramesUI16
SetUIntBuffer, 98	CMcsUsbDacqNet, 191, 192
SetUShortBuffer, 99	ChannelBlock_ReadFramesUI32
SetValue, 99	CMcsUsbDacqNet, 193
WritePipe, 100	ChannelDataEvent
CGilsonDeviceNet, 100	CMcsUsbDacqNet, 219
∼CGilsonDeviceNet, 101	ChannelReset
CGilsonDeviceNet, 101	CMultiBatteryChargerDeviceNet, 299
ConnectSlave, 101	Charge
GetLastAnswer, 101	BatteryState, 30
m_pGilsonDevice, 102	ChargeRegionString
SendBuffered, 101	BatteryState, 30
SendImmediate, 101	ChargeString
SendImmediateGetResponse, 101	BatteryState, 30
CGrapheneFunctionNet, 102	CHiClampDeviceNet, 109
!CGrapheneFunctionNet, 103	CHiClampDeviceNet, 110
·	•
~CGrapheneFunctionNet, 103	RoboDacq, 110
CGrapheneFunctionNet, 103	CHLADacqNet, 110
GetDACOffset, 103, 104	CHLADacqNet, 111
GetVdVs, 104	CHLADeviceNet, 111
GetVdVsDAC, 104, 105	CHLADeviceNet, 111
GetVoltageRange, 105	HLADacq, 112
GetVoltageReached, 105	SerialPort, 112
GetVoltageResolution, 106	CHWInfo
SetDACOffset, 106	CMcsUsbDacqNet::CHWInfo, 113
SetVds, 107	CIntanMea_FunctionNet, 114
SetVdVs, 107	AmplifierSettle, 115
SetVdVsDAC, 108	BeginImpedanceCheck, 115
SetVgs, 108	CIntanMea_FunctionNet, 115
SetVoltageRange, 109	GetDSPHighPassByIndex, 115
SetVoltageResolution, 109	GetImpedanceResult, 115
ChangeSerialNumber	GetIntanRegister, 115
CMcsUsbFactoryNet, 223	GetLowerFrequencyByIndex, 115
channel	GetUpperFrequencyByIndex, 116
CMeaAudioFunctionNet::s_setaudionet, 616	SetBandwidthByIndex, 116
CW2100_FunctionNet::AudioChannelsNet, 29	-
	SetDSPHighPagePulpdox 116
ChannelBlock_AvailFrames	SetDSPHighPassByIndex, 116
CMcsUsbDacqNet, 180, 181	SetIntanRegister, 116
ChannelBlock_ReadAsFrameArrayI16	CInterfaceboard2FunctionNet, 116
CMcsUsbDacqNet, 181	!CInterfaceboard2FunctionNet, 117
ChannelBlock_ReadAsFrameArrayl32	\sim CInterfaceboard2FunctionNet, 117

CInterfaceboard2FunctionNet, 117	IsUserTriggerEnabled, 125
GetloVoltage, 117	ReadClipping, 126
SetloVoltage, 118	SendCommand, 126
CInterfaceboardFunctionNet, 118	SetAdcOffset, 126
!CInterfaceboardFunctionNet, 119	SetAdcOffsetPermanent, 127
\sim CInterfaceboardFunctionNet, 119	SetAudioOutDacParameter, 127
CInterfaceboardFunctionNet, 119	SetDacIdleValue, 127
GetCardinalDacqSamplerate, 119	SetDacOffset, 127
GetCardinalStgOutputrate, 119	SetDacOffsetPermanent, 128
SetCardinalDacqSamplerate, 119	SetDacUseIdleValue, 128
SetCardinalStgOutputrate, 120	SetDigOutState, 128
ClampAmpRestart	SetEEpromPage, 128
CRoboDacqNet, 380	SetSampleInterval, 129
ClearBuffers	StimulusFunction, 129
CMcsUsbDacqNet, 194	CloseAllValves
ClearChannel	CRoboFluidDeviceNet, 406
CDigOutStimulatorFunctionNet, 60	ClosePipe
ClearChannel_PrepareAndSendData	CGenericDevelopDeviceNet, 89
CStg200xDownloadNet, 479	ClosePlateClamp
CStimulusFunctionNet, 486	CMultiwellDeviceNet, 308
ClearChannelData	CMcsBus_AxisParametersNet, 129
CStg200xDownloadBasicNet, 470	~CMcsBus_AxisParametersNet, 130
CStimulusFunctionNet, 486	CMcsBus AxisParametersNet, 129
CW2100_StimulatorFunctionNet, 533	GetAxisParametersSignedEeprom, 130
ClearMultiplexedData	GetAxisParametersUnsignedEeprom, 130
CStimulusFunctionNet, 486	SetAxisParametersEeprom, 130
ClearSTGOutput	CMcsBus_ExtensionNet, 130
CCMOSMea_FunctionNet, 36	~CMcsBus_ExtensionNet, 131
ClearStimulusParametersCache	CMcsBus_ExtensionNet, 131
CW2100_FunctionNet, 526	GetLEDSwitch, 131
ClearSyncData	SetLEDSwitch, 131
CStg200xDownloadBasicNet, 470	CMcsBus_FYIExtensionNet, 131
CStimulusFunctionNet, 486	~CMcsBus_FYIExtensionNet, 132
ClearTable	CMcsBus_FYIExtensionNet, 132
CRetinaLedDeviceNet, 368	GetDIO, 132
ClearTableName	
	GetSingleHeater, 132
CWarnerValveControllerDeviceNet, 561	GetValves, 132
ClearUserDefinedNameCache	SetDIO, 132
CW2100_FunctionNet, 526	SetSingleHeater, 133
ClearValveTable	SetValves, 133
CWarnerValveControllerDeviceNet, 561	CMcsBus_MotorControlNet, 133
CLIH3DeviceNet, 120	~CMcsBus_MotorControlNet, 136
!CLIH3DeviceNet, 122	CMcsBus_MotorControlNet, 136
~CLIH3DeviceNet, 122	GetMCAcceleration, 136
CLIH3DeviceNet, 122	GetMCAccelerationEeprom, 136
DummyCommand, 122	GetMCAccelerationShortCommand, 137
EnableUserTrigger, 122	GetMCAxisRevisionEeprom, 137
ErasePermanentAdcOffset, 123	GetMCBreakCurrent, 137
ErasePermanentDacOffset, 123	GetMCBreakCurrentEeprom, 137
GetAdcOffset, 123	GetMCConfig, 137
GetAudioOutDacParameter, 123	GetMCConfigEeprom, 137
GetDacIdleValue, 124	GetMCCurrent, 137
GetDacOffset, 124	GetMCCurrentEeprom, 138
GetDacqRunStatus, 124	GetMCCurrentMode, 138
GetDacUseIdleValue, 124	GetMCCurrentModeEeprom, 138
GetDigInState, 125	GetMCCurrentModeShortCommand, 138
GetEEpromPage, 125	GetMCCurrentPosition, 138
GetSampleInterval, 125	GetMCCurrentShortCommand, 138

GetMCCurrentSpeed, 138	SetMCReferenceCurrent, 147
GetMCMaxAcceleration, 139	SetMCReferenceCurrentEeprom, 147
GetMCMaxAccelerationEeprom, 139	SetMCRegulatorGain, 148
GetMCMaxCurrent, 139	SetMCRegulatorGainEeprom, 148
GetMCMaxCurrentEeprom, 139	SetMCRotation, 148
	SetMCScalingFactor, 148
GetMCMaxSpeed, 139	G .
GetMCMaxSpeedEeprom, 139	SetMCScalingFactorEeprom, 148
GetMCMaxTravel, 139	SetMCSpeed, 148
GetMCMaxTravelEeprom, 140	SetMCSpeedEeprom, 149
GetMCMaxTravelShortCommand, 140	SetMCSpeedShortCommand, 149
GetMCMovement, 140	SetMCSpeedUnitEeprom, 149
GetMCNewPosition, 140	SetMCStandbyCurrent, 149
GetMCOutputOnOff, 140	SetMCStandbyCurrentEeprom, 149
GetMCPhase, 140	SetMCStandbyTime, 149
GetMCPhaseOffset, 140	SetMCStandbyTimeEeprom, 150
GetMCReference, 141	SetSubChannel, 150
GetMCReferenceCurrent, 141	StartMCMovement, 150
GetMCReferenceCurrentEeprom, 141	StopMCMovement, 150
GetMCRegulatorGain, 141	CMcsBus_SensorNet, 150
GetMCRegulatorGainEeprom, 141	\sim CMcsBus_SensorNet, 152
GetMCScalingFactor, 141	CatchAmpGetAdcMean, 152
GetMCScalingFactorEeprom, 141	CatchAmpGetAdcValue, 153
GetMCSpeed, 142	CatchAmpGetAdcValueH, 153
GetMCSpeedEeprom, 142	CatchAmpGetAdcValueL, 153
GetMCSpeedShortCommand, 142	CatchAmpGetDacAmplitude, 153
GetMCSpeedUnitEeprom, 142	CatchAmpGetDacEnable, 153
GetMCStandbyCurrent, 142	CatchAmpGetDacOffset, 153
GetMCStandbyCurrentEeprom, 142	CatchAmpGetPwmEnable, 153
GetMCStandbyTime, 142	CatchAmpSetDacAmplitude, 153
GetMCStandbyTimeEeprom, 143	CatchAmpSetDacEnable, 154
GetSubChannel, 143	CatchAmpSetDacOffset, 154
SetMCAcceleration, 143	CatchAmpSetPwmEnable, 154
SetMCAccelerationEeprom, 143	CMcsBus_SensorNet, 152
SetMCAccelerationShortCommand, 143	Get2AnalogInput, 154
SetMCAxisRevisionEeprom, 143	Get2DigitalInput, 154
SetMCBreakCurrent, 143	Get4ADC, 154
SetMCBreakCurrentEeprom, 144	Get4ADCAverage, 154
SetMCConfig, 144	Get4ADCCatchampAverageShift, 154
SetMCConfigEeprom, 144	Get4ADCMode, 155
SetMCCurrent, 144	Get4DAC, 155
SetMCCurrentEeprom, 144	GetADCs, 155
SetMCCurrentMode, 144	GetADCsLoop, 155
SetMCCurrentModeEeprom, 145	GetBubbleStatus, 155
SetMCCurrentModeShortCommand, 145	GetDACs, 155
SetMCCurrentPosition, 145	GetDetectionThreshold, 155
SetMCCurrentShortCommand, 145	GetDetectorValue, 155
SetMCMaxAcceleration, 145	GetLatency, 156
SetMCMaxAccelerationEeprom, 145	GetLatencyCounter, 156
SetMCMaxCurrent, 146	GetMinimalThreshold, 156
SetMCMaxCurrentEeprom, 146	GetMovePump, 156
SetMCMaxSpeed, 146	GetPiezoState, 156
SetMCMaxSpeedEeprom, 146	GetPressure, 156
SetMCMaxTravel, 146	GetPressureOffset, 157
SetMCMaxTravelEeprom, 146	GetRegulationTimeouts, 157
SetMCMaxTravelShortCommand, 147	GetRegulatorFactor, 157
SetMCNewPosition, 147	GetRegulatorOnOff, 157
SetMCOutputOnOff, 147	GetRegulatorStatus, 157
SetMCReference, 147	GetRotatePump, 157
OCTIVIOLICICICIOS, 147	deli iolalei ullip, 137

0.10 1.5 1.1 1.5	
GetSamplePeriode, 158	CMcsBusNet::GetMode, 168
GetSollPressure, 158	CMcsBusNet::GetModeEeprom, 168
GetSyncState, 158	CMcsBusNet::SetMode, 169
Set4ADCCatchampAverageShift, 158	CMcsBusNet::SetModeEeprom, 169
Set4ADCMode, 158	GetBusAddress, 169
Set4DAC, 158	GetBusAddressEeprom, 169
SetDACs, 158	GetCommand, 169, 170
SetDetectionThreshold, 159	GetHWRevisionEeprom, 170
SetLatency, 159	SetBusAddress, 170
SetMinimalThreshold, 159	SetBusAddressEeprom, 170
SetMovePump, 159	SetCommand, 170, 171
SetPiezoState, 159	SetHWRevisionEeprom, 171
SetPressureOffset, 159	CMcsBusNet::GetMode
SetRegulationTimeouts, 159	CMcsBusNet, 168
SetRegulatorFactor, 160	CMcsBusNet::GetModeEeprom
SetRegulatorOnOff, 160	CMcsBusNet, 168
SetRotatePump, 160	CMcsBusNet::SetMode
SetSamplePeriode, 160	CMcsBusNet, 169
SetSollPressure, 160	CMcsBusNet::SetModeEeprom
StartSync, 160	CMcsBusNet, 169
TactSwitchGetState, 160	
,	CMcsUsbDacqNet, 171
TactSwitchSetDisplay, 161	~CMcsUsbDacqNet, 177
CMcsBus_TempSensorNet, 161	AddSelectedChannelsQueue, 177, 179, 180
~CMcsBus_TempSensorNet, 161	ChannelBlock_AvailFrames, 180, 181
CMcsBus_TempSensorNet, 161	ChannelBlock_ReadAsFrameArrayl16, 181
GetNanoVoltsPerKelvin, 162	ChannelBlock_ReadAsFrameArrayl32, 182, 183
GetTemperatur, 162	ChannelBlock_ReadAsFrameArrayUI16, 183, 184
GetThermoOffset, 162	ChannelBlock_ReadAsFrameArrayUl32, 185, 186
GetThermoTemp, 162	ChannelBlock_ReadFramesDictl16, 186
GetThermoVoltage, 162	ChannelBlock_ReadFramesDictl32, 187
SetNanoVoltsPerKelvin, 162	ChannelBlock_ReadFramesDictUI16, 187
SetThermoOffset, 162	ChannelBlock_ReadFramesDictUl32, 188
CMcsBus_VoltageModeNet, 163	ChannelBlock_ReadFramesI16, 189
\sim CMcsBus_VoltageModeNet, 164	ChannelBlock_ReadFramesl32, 190, 191
CMcsBus_VoltageModeNet, 164	ChannelBlock_ReadFramesUI16, 191, 192
GetVMMaxNegativeCurrent, 164	ChannelBlock_ReadFramesUl32, 193
GetVMMaxNegativeCurrentEeprom, 164	ChannelDataEvent, 219
GetVMMaxNegativeVoltage, 164	ClearBuffers, 194
GetVMMaxNegativeVoltageEeprom, 164	CMcsUsbDacqNet, 177
GetVMMaxPositiveCurrent, 165	CMcsUsbDacqNet::GetFilterProperties, 194
GetVMMaxPositiveCurrentEeprom, 165	Error_Callback_Aquisition_Stopped, 218
GetVMMaxPositiveVoltage, 165	Error_Callback_Data_lost, 218
GetVMMaxPositiveVoltageEeprom, 165	Error_Callback_Frames_Lost, 218
GetVMOutputOnOff, 165	Error_Callback_Packet_Error, 218
GetVMVoltage, 165	Error_Callback_Queue_Full, 219
SetVMMaxNegativeCurrent, 165	Error_Callback_RingQueue_Full, 219
SetVMMaxNegativeCurrentEeprom, 166	ErrorEvent, 219
SetVMMaxNegativeVoltage, 166	GetAdapterType, 194
SetVMMaxNegativeVoltageEeprom, 166	GetAdcDataFormat, 194
SetVMMaxPositiveCurrent, 166	GetAdcZero, 194
SetVMMaxPositiveCurrentEeprom, 166	GetAnalogValueUnit, 195
SetVMMaxPositiveVoltage, 166	GetChannelDataFillSize, 195
SetVMMaxPositiveVoltageEeprom, 167	GetChannelLayout, 195
SetVMOutputOnOff, 167	GetChannelsInBlock, 195
SetVMVoltage, 167	GetDataFormat, 195
CMcsBusNet, 167	GetDataMode, 195
~CMcsBusNet, 168	GetDataMode, 193 GetDigitalSource, 196, 197
CMcsBusNet, 168	GetFilterProperty, 198
OMOSDUSINGI, 100	deli illeri roperty, 130

GetGroupChannelDatal16, 198	CMcsUsbDeviceStatePushNet, 221
GetGroupChannelDatal32, 198	McsUsbDeviceStateEvent, 221
GetGroupChannelDataUI16, 199	TriggerStatus, 221
GetGroupChannelDataUl32, 199	CMcsUsbFactoryNet, 221
GetHardwareMaxRange, 200	~CMcsUsbFactoryNet, 223
GetHardwareMinRange, 200	ChangeSerialNumber, 223
GetMaxSamplingFrequency, 200	CMcsUsbFactoryNet, 223
GetMeaLayout, 200	Coldstart, 223
GetMinSamplingFrequencyStepsize, 201	CompareFirmware, 224
GetNumberOfDataBits, 201	DownloadFirmware, 224
GetPoti, 201	FindFirmwareVersionMagicInBuffer, 224
GetResolutionPerDigit, 201	FX3MCSDataAddress, 229
GetSamplerate, 201	FX3MCSDataDeviceIdOffset, 229
GetVoltageRangeIndex, 201	FX3MCSDatalFB1ImageOffset, 229
GetVoltageRangeInMicroVolt, 202	FX3MCSDatalFB2ImageOffset, 229
GetVoltageRangeInMilliVolt, 202	FX3MCSDataVersionOffset, 229
HWInfo, 202	GetChecksumFromFX3Image, 224
Samplerate, 219	GetDestination, 224
SendStartDacq, 202	GetDestinationDisplayLabel, 224
SendStartStgAndDacq, 202	GetDestinationName, 224, 225
SendStopDacq, 203	GetDestinationName, 224, 225 GetDestinationSerialNumber, 225
SendStopStgAndDacq, 203	
	GetDestinationTargetAddress, 225
SendStopStgAndDacqWithOptions, 203	GetFirmwareVersionFromFile, 225
SetDataMode, 204	GetFirmwareVersionFromHexFile, 225
SetDigitalSource, 204–206	GetNumDestinations, 225
SetPoti, 206	GetUSBDeviceIDFromFX3Image, 226
SetSamplerate, 206	GetUsercodeFromBitFile, 226
SetSelectedChannels, 207–209	GetUsercodeFromFlash, 226
SetSelectedChannelsQueue, 209–211	GetXilinxFlashOffset, 226
SetSelectedData, 211–213	GetXilinxFlashReadCommand, 226
SetupGroupDacqQueue, 213	LoadUserFirmware, 226, 227
SetVoltageRangeByIndex, 213	ReadBlockFromFlash, 227
SetVoltageRangeInMicroVolt, 214	ReadBlockFromIFBGlobalEEprom, 227
StartDacq, 214, 215	ReadBlockFromNVMEM, 227
StartLoop, 216, 217	SetDestinationSerialNumber, 227
StopDacq, 218	UpdateFirmware, 227–229
StopLoop, 218	CMcsUsbFunctionNet, 230
CMcsUsbDacqNet::CHWInfo, 112	!CMcsUsbFunctionNet, 231
CHWInfo, 113	\sim CMcsUsbFunctionNet, 231
GetAvailableSampleRates, 113	CMcsUsbFunctionNet, 230, 231
GetAvailableVoltageRangesInMicroVolt, 113	m_pMcsUsb, 231
GetAvailableVoltageRangesInMicroVoltAnd-	m_pMcsUsbFunction, 231
StringsInMilliVolt, 113	ThrowCUsbExceptionNetOnError, 231
GetNumberOfHWADCChannels, 113	CMcsUsbFunctionPointerContainer, 231
GetNumberOfHWDigitalChannels, 113	CMcsUsbListEntryNet, 231
IsDigitalChannelDedicated, 114	~CMcsUsbListEntryNet, 232
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet,	Deviceld, 236
523	DeviceName, 236
CVoltageRangeInfoNet, 523	Equals, 233
VoltageRangeDisplayStringMilliVolt, 523	GetEntry, 234
VoltageRangeInMicroVolt, 523	GetEntryCount, 235
CMcsUsbDacqNet::GetFilterProperties	HwVersion, 236
CMcsUsbDacqNet, 194	Manufacturer, 236
CMcsUsbDeviceStatePushFunctionNet, 219	Product, 236
CMcsUsbDeviceStatePushFunctionNet, 220	SerialNumber, 236
McsUsbDeviceStateEvent, 220	SetStringFormat, 235
TriggerStatus, 220	ToString, 236
CMcsUsbDeviceStatePushNet, 220	CMcsUsbListNet, 237

!CMcsUsbListNet, 238	ReadRegister, 253
\sim CMcsUsbListNet, 238	ReadRegister32, 253
CMcsUsbListNet, 237	ReadRegisterTimeSlot, 253
Count, 239	RemoveSoftwareKey, 254
DeviceArrival, 239	RescanHeadstage, 254
DeviceRemoval, 239	SerialNumber, 261
GetNumberOfDevices, 238	SetConfiguration, 254
GetUsbListEntries, 238	SetSoftwareKey, 254
GetUsbListEntry, 238	Status_AlreadyConfigured, 256
IsDeviceTypeOf, 239	Status_BadStartFrame, 257
SetStringFormat, 239	Status_Btstuff, 257
CMcsUsbNet, 240	Status_BufferOverrun, 257
!CMcsUsbNet, 244	Status_BufferUnderrun, 257
∼CMcsUsbNet, 244	Status_Canceled, 257
AddSoftwareKey, 244	Status_Canceling, 257
AssociateToThis, 244	Status_ConnectedPipes, 257
CMcsUsbNet, 244	Status_ControlNotOwned, 257
Connect, 244, 245	Status Crc, 257
Disconnect, 246	Status_DataOverrun, 257
EmptyKey, 246	Status_DataToggleMismatch, 257
EnableExceptions, 246	Status DataUnderrun, 258
EraseEepromRegisterPreconfig, 246	Status_DeviceLocked, 258
GetConfiguration, 247	Status_DeviceNotFound, 258
GetDeviceCannotStallOutRequests, 247	Status DeviceRemoved, 258
GetDeviceCapableSpeed, 247	Status_DevNotResponding, 258
GetDeviceOapableOpeed, 247	Status_EndpointHalted, 258
GetDeviceLitatii, 247	Status_ErrorBusy, 258
GetDevicerd, 247 GetDeviceRootHubVendorEnum, 247	Status_ErrorShortTransfer, 258
GetDeviceRootHubVendorID, 247	Status Fifo, 258
GetDeviceRootHubVendorName, 247	Status FrameControlOwned, 258
GetDeviceRootHubVerldofName, 247 GetDeviceSpeed, 248	-
•	Status_InternalHcError, 258 Status InvalidParameter, 259
GetErrorText, 248 GetFirmwareVersion, 248	-
•	Status_Invalid light-unstion 259
GetHardwareRevision, 248	Status_InvalidUrbFunction, 259
GetHeadstageActive, 249	Status_loPending, 259
GetHeadstageID, 249	Status_loTimeout, 259
GetHeadstagePresent, 249	Status_IsochRequestFailed, 259
GetIdent, 249	Status_LastUsbErrorMismatch, 259
GetLastUSBError, 250	Status_NoBandwidth, 259
GetMea21UsbPort, 250	Status_NoMemory, 259
GetNumConfigurations, 250	Status_NoSuchDevice, 259
GetSerialNumber, 250	Status_NotAccessed, 259
GetSoftwareKey, 250	Status_NotSupported, 260
GetSoftwareKeyString, 250	Status_PidCheckFailure, 260
GetStatus, 250	Status_PipeNotLinked, 260
GetStatusOfLastCommand, 251	Status_RequestFailed, 260
GetUsbListEntry, 251	Status_RequestMutexFailed, 260
GetVersion, 251	Status_RequestMutexTimeout, 260
HasSoftwareKey, 251	Status_Stall, 260
IsConnected, 251	Status_Unconfigured, 260
IsDeviceHighSpeed, 251	Status_UnexpectedPid, 260
IsDeviceHighSpeedCapable, 252	ThrowCUsbExceptionNetOnError, 254
IsExceptionsEnabled, 252	TxnGetSerialNumber, 254
MultibootGetCypressImageId, 252	TxnSetSerialNumber, 254
MultibootGetImageId, 252	TxnTestMemoryReadAndCheck, 254
MultibootGetSelectedImage, 252	TxnTestMemoryWrite, 255
MultibootSelectImage, 252	ValidKey, 255
ReadEepromRegisterPreconfig, 253	WPAError_ScanningIsPending, 260

WriteEepromRegisterPreconfig, 255	SetSlope, 275
WriteRegister, 255, 256	Start, 275
WriteRegister32, 256	Stop, 276
WriteRegisterArray, 256	CMeaDeviceNet, 276
WriteRegisterTimeSlot, 256	~CMeaDeviceNet, 278
WriteRegisterValue, 256	AnalogGain, 283
CMcsUsbPointerContainer, 261	CMeaDeviceNet, 277
CMEA2100x256FunctionNet, 261	EnableChecksum, 278
!CMEA2100x256FunctionNet, 262	EnableDigitalIn, 278, 279
~CMEA2100x256FunctionNet, 262	EnableTimestamp, 279
CMEA2100x256FunctionNet, 261, 262	Gain, 283
GetLayoutConfiguration, 262	GetAnalogGain, 280
SetLayoutConfiguration, 262	GetEnumerationSpeed, 280
CMeaAudioFunctionNet, 263	GetGain, 280
CMeaAudioFunctionNet, 263	MeaAudioFunctionNet, 283
GetAudioChannels, 264, 265	MeaDigitalDataFunctionNet, 283
GetNumberOfAudioChannels, 265	MeaFeedbackFunctionNet, 283
SetAudioChannels, 265, 266	MeFunctionNet, 284
CMeaAudioFunctionNet::s setaudionet, 616	SetDigitalOut, 280
amplification, 616	SetNumberOfAnalogChannels, 280
channel, 616	SetNumberOfChannels, 281, 282
CMeaCleanDeviceNet, 266	SetTriggerMaskValue, 282
!CMeaCleanDeviceNet, 268	SetTriggerPeriod, 283
~CMeaCleanDeviceNet, 267	W2100 FunctionNet, 284
CMeaCleanDeviceNet, 267	WClassicFunctionNet, 284
GetCycle, 268	CMeaDigitalDataFunctionNet, 284
GetCycles, 268	CMeaDigitalDataFunctionNet, 284
GetMaxVoltage, 268	GetDigitalData, 285
GetMinVoltage, 268	SetDigitalData, 285
GetOutputVoltage, 268	CMeaFeedbackFunctionNet, 286
GetOutputVoltage, 268	CMeaFeedbackFunctionNet, 286 CMeaFeedbackFunctionNet, 287
	CMeaFeedbackFunctionNet, 287
GetOutputVoltage, 268 GetSlope, 269	
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 CMeaCoatDeviceNet, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetMaxCurrent, 272 GetMaxCurrent, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfRateDetectors, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetMaxCurrent, 272 GetMaxCurrent, 272 GetOffsetCurrent, 273	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIIRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfRateDetectors, 289 FeedbackSetNumberOfSpikeDetectors, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272 GetOffsetCurrent, 273 GetPauseDuration, 273	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetUlRFilterParameter, 288 FeedbackSetNumberOfcagics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCurrentCycle, 272 GetDuration, 272 GetOtputCurrent, 273 GetPauseDuration, 273 GetSlope, 273	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272 GetMaxCurrent, 273 GetOutputCurrent, 273 GetPauseDuration, 273 GetSlope, 273 GetTimeInPause, 273	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetUgic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289 FeedbackSetSpikeDetectorThreshold, 289
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272 GetMaxCurrent, 273 GetOutputCurrent, 273 GetPauseDuration, 273 GetSlope, 273 GetTimeInPause, 273 GetTimeInPlateau, 274	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetUgic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289 FeedbackSetSpikeDetectorThreshold, 289 FeedbackSetTrigger, 290
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272 GetMaxCurrent, 273 GetOutputCurrent, 273 GetPauseDuration, 273 GetTimeInPause, 273 GetTimeInPlateau, 274 IsRunning, 274	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetUgic, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289 FeedbackSetSpikeDetectorThreshold, 289 FeedbackSetTrigger, 290 CMealmpedanceDeviceNet, 290
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCycles, 272 GetDuration, 272 GetMaxCurrent, 273 GetOutputCurrent, 273 GetPauseDuration, 273 GetTimeInPause, 273 GetTimeInPlateau, 274 IsRunning, 274 SetCycles, 274	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter32, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetIlRFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetMkFilter, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfRateDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289 FeedbackSetSpikeDetectorThreshold, 289 FeedbackSetTrigger, 290 CMeaImpedanceDeviceNet, 290 ~CMeaImpedanceDeviceNet, 290
GetOutputVoltage, 268 GetSlope, 269 IsRunning, 269 SetCycles, 269 SetMaxVoltage, 269 SetMinVoltage, 270 SetSlope, 270 Start, 270 Stop, 270 CMeaCoatDeviceNet, 270 !CMeaCoatDeviceNet, 272 ~CMeaCoatDeviceNet, 272 GetCurrentCycle, 272 GetCurrentCycle, 272 GetDuration, 272 GetMaxCurrent, 272 GetOtfsetCurrent, 273 GetPauseDuration, 273 GetSlope, 273 GetTimeInPause, 273 GetTimeInPlateau, 274 IsRunning, 274 SetCycles, 274 SetDuration, 274	CMeaFeedbackFunctionNet, 287 FeedbackGetSampleTimerCount, 287 FeedbackSetAnalogSource, 287 FeedbackSetChannelFilter, 287 FeedbackSetDigitalMapping, 287 FeedbackSetFeedback, 287 FeedbackSetFilterOff, 288 FeedbackSetFilterParameter, 288 FeedbackSetFilterParameter, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetGlobalChannelFilter, 288 FeedbackSetURFilterParameter, 288 FeedbackSetLogic, 288 FeedbackSetNumberOfLogics, 289 FeedbackSetNumberOfRateCounter, 289 FeedbackSetNumberOfRateDetectors, 289 FeedbackSetNumberOfSpikeDetectors, 289 FeedbackSetNumberOfTriggers, 289 FeedbackSetRateCounter, 289 FeedbackSetRateCounter, 289 FeedbackSetRateDetector, 289 FeedbackSetRateDetectorThreshold, 289 FeedbackSetTrigger, 290 CMeaImpedanceDeviceNet, 290 CMeaImpedanceDeviceNet, 290 CMeaImpedanceDeviceNet, 290

GetReady, 291	\sim CMultiwellDeviceNet, 307
GetResult, 291	ClosePlateClamp, 308
SetImpedanceTestFrequency, 291	CMultiwellDeviceNet, 307
StartMeasurement, 291	GetPlateClampLockState, 308
CMeasureTableDeviceNet, 291	GetPlateClampState, 308
CMeasureTableDeviceNet, 292	GetPlateMux, 308, 309
Sensor, 292	GetPlateType, 309
CMeaSwitchDeviceNet, 292	GetPowerMuxPlate, 309
~CMeaSwitchDeviceNet, 293	IsPlateTypeValid, 310
CMeaSwitchDeviceNet, 293	LockPlateClamp, 310
GetNumber, 293	OpenPlateClamp, 310
GetPattern, 293	SetPlateMux, 310, 311
GetPatternBool, 294	SetPlateType, 311
SetPattern, 294	SetPowerMuxPlate, 311
SetPatternBool, 294	StopPlateClamp, 312
CMeaUSBDeviceNet, 294	UnlockPlateClamp, 312
~CMeaUSBDeviceNet, 295	CMultiwellOptoStimFunctionNet, 312
CMeaUSBDeviceNet, 295	!CMultiwellOptoStimFunctionNet, 313
CMeFunctionNet, 295	~CMultiwellOptoStimFunctionNet, 313
!CMeFunctionNet, 296	CMultiwellOptoStimFunctionNet, 313
~CMeFunctionNet, 296	GetAbsMaxCurrentInMicroAmp, 313
CMeFunctionNet, 296	GetColorRgb, 313
SetTransformer, 296	GetColorStr, 314
CMosMea	GetMaxDurationHighCurrentInMicroSec, 314
CCMOSMeaDeviceNet, 47	GetMaxDutyCycleHighCurrent, 314
CMultiBatteryChargerDeviceNet, 297	GetPermanentCurrentInMicroAmp, 315
!CMultiBatteryChargerDeviceNet, 298	GetWaveLengthInNanometer, 315
~CMultiBatteryChargerDeviceNet, 298	SetAbsMaxCurrentInMicroAmp, 315
CapacityTest, 298	SetColorRgb, 315
ChannelReset, 299	SetColorStr, 316
CMultiBatteryChargerDeviceNet, 298	SetMaxDurationHighCurrentInMicroSec, 316
GetBatteryVoltage, 299	SetMaxDutyCycleHighCurrent, 316
GetChannels, 299	SetPermanentCurrentInMicroAmp, 316
GetChannelState, 299	SetWaveLengthInNanometer, 317
GetChargeCapacity, 300	CNF_GenDeviceNet, 317
GetChargeCurrent, 300	~CNF_GenDeviceNet, 317
GetChargingMode, 300	CNF_GenDeviceNet, 317
GetChargingPCoefficient, 300	Set_Values, 317
GetDischargeCapacity, 301	COctoPotDeviceNet, 318
GetDischargeCurrent, 301	BurnAdcOffset, 319
GetDischargeCurrentSetPoint, 301	BurnDacOffset, 319
GetFinalDischargeVoltage, 302	COctoPotDeviceNet, 318, 319
GetRatedCapacity, 302	EnableChecksum, 319
SetChargingMode, 302	EnableDigitalIn, 319
SetChargingPCoefficient, 302	EnableTimestamp, 319
SetDischargeCurrentSetPoint, 303	GetAdcOffset, 319
SetFinalDischargeVoltage, 303	GetDacOffset, 319
SetRatedCapacity, 303	PatternListStart, 319
SetRatedCapacityVolatile, 303	RampStart, 320
CMultiwellCallbackFunctionNet, 304	ResetAdcOffset, 320
!CMultiwellCallbackFunctionNet, 305	ResetDacOffset, 320
~CMultiwellCallbackFunctionNet, 305	SetAdcOffset, 320
CMultiwellCallbackFunctionNet, 305	SetAmplificationSwitch, 320
GetPlateClampStateByHeadstage, 305	SetBathclamp, 320
GetPlateClampStateByHeadstageEvent, 306	SetChannelSwitch, 320
OnGetPlateClampStateByHeadstage, 305	SetDacAutoControl, 320
CMultiwellDeviceNet, 306	SetDacOffset, 320
!CMultiwellDeviceNet, 307	SetDacValue, 321

SetNumberOfChannels, 321	CPatchServerDeviceNet, 326
SetOutputRate, 321	Sensor, 326
SetPatternListEntry, 321	CPathIdentDeviceNet, 326
SetPidParameter, 321	\sim CPathIdentDeviceNet, 327
SetRampParameter, 321	CPathIdentDeviceNet, 327
SetSineParameter, 321	Measure, 327
SineStart, 322	Set_Values, 327
COkuvisionStimulatorDeviceNet, 322	CPedoterDeviceNet, 327
~COkuvisionStimulatorDeviceNet, 323	!CPedoterDeviceNet, 328
COkuvisionStimulatorDeviceNet, 323	~CPedoterDeviceNet, 328
GetCheckVoltage, 323	CPedoterDeviceNet, 328
GetCurrentFactor, 323	GetCommand, 328
GetDACOffset, 323	SetCommand, 328
GetMaxPower, 323	CPeristalticPumpDeviceNet, 329
GetMaxVoltage, 323	\sim CPeristalticPumpDeviceNet, 329
GetPulseform, 323	CPeristalticPumpDeviceNet, 329
GetRTC, 324	McsBus MotorControl, 330
GetStimulatorStatus, 324	CPgaDeviceNet, 330
GetVoltage, 324	~CPgaDeviceNet, 330
SetCheckVoltage, 324	ApplyGains, 331
-	• • •
SetCurrentFactor, 324	CPgaDeviceNet, 330
SetDACOffset, 324	DefineAmplification, 331
SetMaxPower, 325	DefineFrequencyRange, 331
SetMaxVoltage, 325	DefineNumAmplifications, 331
SetPulseform, 325	DefineNumFrequencyRanges, 331
SetRTC, 325	GetAmplification, 331
Coldstart	GetFrequencyRange, 331
CMcsUsbFactoryNet, 223	GetGain, 332
CommaPositionA	GetNumAmplifications, 332
CFilterCoefficientsNet::s_FilterAttributesNet, 616	GetNumFrequencyRanges, 332
CommaPositionB	SetGain, 332
CFilterCoefficientsNet::s_FilterAttributesNet, 616	CPositionIIDeviceNet, 332
CompareFirmware	!CPositionIIDeviceNet, 334
CMcsUsbFactoryNet, 224	\sim CPositionIIDeviceNet, 334
CompareTo	CPositionIIDeviceNet, 334
HeadStageIDType, 601	GetCoilCommunication, 334
CompensateElectrodeOffset	GetDebugData, 334
CWarnerUssingFunctionNet, 543	GetEventData, 335
Connect	GetImplantCurrentSetpoint, 335
CMcsUsbNet, 244, 245	GetImplantResult, 336
CRFFunctionNet, 370	GetImplantState, 336
ConnectDevice	GetOnOff, 336
CRadioControledDevicesNet, 365	GetPowerStrength, 337
	GetRTC, 337
ConnectedImp	
CPositionImpDeviceNet, 342	GetStateDebugData, 337
ConnectImp	GetStateEventData, 338
CPositionImpDeviceNet, 342	RFFunction, 341
ConnectSlave	SetImplantCurrentSetpoint, 338
CGilsonDeviceNet, 101	SetPowerStrength, 338
ControlState	SetRTC, 338
HeadStageIDTypeState, 604	SetStateDebugData, 340
CornerFrequency	SetStateEventData, 340
CFilterPropertyNet, 72	SwitchOnOff, 340
CornerFrequencymHz	CPositionImpDeviceNet, 341
CFilterPropertyNet, 73	!CPositionImpDeviceNet, 342
	·
Count CM and Jahl Jahlah 2000	~CPositionImpDeviceNet, 342
CMcsUsbListNet, 239	ConnectedImp, 342
CPatchServerDeviceNet, 325	ConnectImp, 342

	CPositionImpDeviceNet, 342	GetUseBubble, 357
	GetDeviceList, 342	SetAnalogVoltages, 357
	GetImpld, 343	SetPumpCouple, 358
	GetRFFrequency, 343	SetPumpEnableSpeedRatio, 358
	SetDeviceList, 343	SetPumpFastOnOff, 358
	SetImpld, 343	SetPumpFastSpeed, 358
	SetRFFrequency, 344	SetPumpFunctionSpeeds, 358
CPF	CDeviceNet, 344	SetPumpManualOnOff, 358
	CPPCDeviceNet, 344	SetPumpMaxSpeed, 358
	McsBus, 345	SetPumpModeType, 358
	McsBus_MotorControl, 345	SetPumpSpeedRatio, 359
	McsBus_Sensor, 345	SetPumpSpeedUnit, 359
	PPCFunction, 345	SetUseBubble, 359
CPF	CFunctionNet, 345	CPPSDeviceNet, 359
	!CPPCFunctionNet, 347	CPPSDeviceNet, 360
	~CPPCFunctionNet, 347	CProgramPressureCurveNet, 360
	CPPCFunctionNet, 346	!CProgramPressureCurveNet, 361
	FirePressurePulse, 347	~CProgramPressureCurveNet, 360
	GetAnalogVoltage, 347	CProgramPressureCurveNet, 360
	GetAnalogVoltageRange, 347	GetRepeats, 361
	GetDigitalIn, 349	Program, 361
	GetPressureRange, 349	SetRepeats, 361
	GetPumpModeType, 349	CPulseGeneratorFunctionNet, 361
	GetPumpSpeedUnit, 350	!CPulseGeneratorFunctionNet, 362
	GetSupplyVoltage, 350	~CPulseGeneratorFunctionNet, 362
	GetValveActive, 350	CPulseGeneratorFunctionNet, 362
	IsBusy, 350	GetModeSelect, 363
	LoadPressure, 352	GetPeriod, 363
	MeasureReservoir, 352	GetPulseLength, 363
	SetAnalogVoltageRange, 352	SetModeSelect, 364
	SetPressureOffset, 352	SetPeriod, 364
	SetPressureRange, 352	SetPulseLength, 364
	SetPumpModeType, 353	CRadioControledDevicesNet, 364
	SetPumpSpeedUnit, 353	ConnectDevice, 365
	SetValveActive, 353	CRadioControledDevicesNet, 365
CPF	PS_DeviceNet, 354	DisConnectDevice, 365
_	CPPS DeviceNet, 354	GetDeviceNames, 365
	McsBus, 354	GetFrequency, 366
	McsBus_MotorControl, 354	HasRadioControl, 366
	McsBus Sensor, 354	SetFrequency, 366
	PPS_Function, 354	StillConnected, 366
CPF	S_FunctionNet, 355	CreateSideband
	CPPS FunctionNet, 355, 356	CStimulusFunctionNet, 486
	GetAnalogVoltage, 356	CreateWirelessHeadstageSerialNumberString
	GetAnalogVoltages, 356	CWirelessBaseFunctionNet, 588
	GetBubbleState, 356	CRegionOfInterestRect
	GetDigitalIn, 356	CCMOSMeaDeviceNet::CRegionOfInterestRect,
	GetPumpCouple, 356	366
	GetPumpEnableSpeedRatio, 356	CRetinaLedDeviceNet, 367
	GetPumpFastOnOff, 356	~CRetinaLedDeviceNet, 368
	GetPumpFastSpeed, 356	AddLoopEntry, 368
	GetPumpFunctionSpeeds, 357	AddTableEntry, 368
	GetPumpManualOnOff, 357	ClearTable, 368
	GetPumpMaxSpeed, 357	CRetinaLedDeviceNet, 368
	GetPumpModeType, 357	GetTablepointer, 368
	GetPumpSpeedRatio, 357	SetLED, 368
	GetPumpSpeedUnit, 357	SetLumi, 368
	GetSupplyVoltage, 357	SetPersistency, 368
	The Arms of Arms	J ;

SetRepeat, 369	GetAllDigout, 381
SetTablepointer, 369	GetCapacityC, 381
SetTrigger, 369	GetCapacityV, 381
CRFFunctionNet, 369	GetCapacityX, 382
!CRFFunctionNet, 370	GetClampAmpSerialNumber, 382
~CRFFunctionNet, 370	GetCommand, 382
Connect, 370	GetConfigurationBit, 382
CRFFunctionNet, 370	GetConfigurationBitAxc, 382
GetAvailableDeviceList, 371	GetConfigurationBitBlu_Led, 382
GetAvailableDeviceListEx, 371	GetConfigurationBitBlu LedToggleFast, 382
GetAvailableStateList, 371	GetConfigurationBitBlu LedToggleSlow, 382
GetAvailableStateListEx, 371	GetConfigurationBitCC Gen, 382
GetBaseFrequency, 372	GetConfigurationBitCV_Gen, 382
GetConnectedDevice, 372	GetConfigurationBitRC_Gen, 382
GetState, 372	GetConfigurationBitRed_Led, 383
	-
GetTestMode, 372	GetConfigurationBitRed_LedSaturation, 383
GetWorkingFrequency, 373	GetConfigurationBitRed_LedToggleFast, 383
SetBaseFrequency, 373	GetConfigurationBitRed_LedToggleSlow, 383
SetTestMode, 373	GetConfigurationBitRelais, 383
SetWorkingFrequency, 373	GetConfigurationBitRV_Gen, 383
CRobo_FYIProgram_FunctionNet, 374	GetConfigurationBits, 383
CRobo_FYIProgram_FunctionNet, 374	GetConfigurationBitStream, 383
GetLength, 374	GetConfigurationBitSupply, 383
GetState, 375	GetCrossTalkOffset, 383
GetValve1, 375	GetCrossTalkOptimum, 383
GetValve2, 375	GetDigout, 384
SetLength, 375	GetDisplayText, 384
SetValve1, 375	GetDownsampleFactor, 384
SetValve2, 375	GetFilter, 384
Start, 375	GetFilterCoeffs, 384
CRobo_FYITemp_FunctionNet, 375	GetIC, 384
CRobo_FYITemp_FunctionNet, 376	GetlClamp, 384
GetlCoeff, 376	GetICOffset, 384
GetMaxPower, 376	GetlGain, 384
GetPCoeff, 376	GetNIC_MS, 384
GetRegulatorOnOff, 376	GetNUC_MS, 384
GetSollTemp, 376	GetNUV_MS, 385
SetlCoeff, 377	GetPGain, 385
SetMaxPower, 377	GetRecordingNumber, 385
SetPCoeff, 377	GetResistanceC, 385
SetRegulatorOnOff, 377	GetResistanceV, 385
SetSollTemp, 377	GetScreen, 385
CRoboDacqNet, 377	GetSimulation, 385
CancelTableLoop, 380	GetUC, 385
CancelTableLoopAndStopTable, 380	GetUClamp, 385
ClampAmpRestart, 380	GetUCOffset, 385
CRoboDacqNet, 380	GetUpdateDisplay, 385
DoRamp, 380	GetUV, 386
Emu GetCellCapacity, 380	GetUVOffset, 386
Emu_GetCellPotential, 380	GetXGain, 386
Emu_GetCellResists, 380	RunTable, 386
Emu_GetElectrodeResists, 381	SetAllDigout, 386
	-
Emu_GetNoise, 381	SetCommand, 386
Emu_SetCellCapacity, 381	SetConfigurationBit Ava 386
Emu_SetCellPotential, 381	SetConfigurationBitAxc, 386
Emu_SetCellResists, 381	SetConfigurationBitBlu_Led, 386
Emu_SetBleicrodeResists, 381	SetConfigurationBitBlu_LedToggleFast, 387
Emu_SetNoise, 381	SetConfigurationBitBlu_LedToggleSlow, 387

0 - + 0 + i Di+ 0 0	0.07	0-+5
SetConfigurationBitCC		GetErrorVoltage5V, 396
SetConfigurationBitCV		GetErrorVoltageAirvalve, 396
SetConfigurationBitRC		GetErrorVoltageRs485A, 396
SetConfigurationBitRe		GetErrorVoltageRs485B, 396
SetConfigurationBitRe		GetErrorVoltageValves, 396
SetConfigurationBitRe	d_LedToggleFast, 387	GetInMovement, 396
SetConfigurationBitRe	d_LedToggleSlow, 387	GetMinPressure, 396
SetConfigurationBitRe	lais, 387	GetMovementError, 396
SetConfigurationBitRV	⁷ Gen, 388	GetVoltage12V, 397
SetConfigurationBitStr		GetVoltage12VLimit, 397
SetConfigurationBitSu		GetVoltage5V, 397
SetCrossTalkOffset, 38	• • •	GetVoltage5VLimit, 397
SetCrossTalkOptimum		GetVoltageAirvalve, 397
SetDigout, 388	, 000	GetVoltageAirvalveLimit, 397
SetDisplayText, 388		GetVoltageRs485A, 397
• •	× 000	-
SetDownsampleFactor	1, 300	GetVoltageRs485ALimit, 397
SetFilter, 388		GetVoltageRs485B, 397
SetFilterCoeffs, 389		GetVoltageRs485BLimit, 398
SetlClamp, 389		GetVoltageValves, 398
SetICOffset, 389		GetVoltageValvesLimit, 398
SetlGain, 389		McsBus, 405
SetNoFilterCoeffs, 389	€	McsBus_MotorControl, 405
SetPGain, 389		McsBus_XY, 402
SetRecordingNumber,	, 389	McsBus_ZI, 402
SetScreen, 389		MoveAbs, 398
SetSimulation, 389		RoboError_AnotherMaster, 402
SetUClamp, 390		RoboError_Base, 402
SetUCOffset, 390		RoboError_CannotEscapeEndSwitch, 402
SetUVOffset, 390		RoboError_CommandAlreadyInProgress, 402
SetXGain, 390		RoboError_CommandNotPossible, 402
StopTable, 390		RoboError_CommunicationTimeout, 402
		RoboError_DacqNotReady, 402
Table_Wait, 390		
TableDefBegin, 390		RoboError_DLLMovementTimeout, 403
TableDefEnd, 390		RoboError_FindReferenceMethod, 403
UpdateDisplay, 390		RoboError_GilsonCommandPending, 403
CRoboDeviceNet, 391		RoboError_GilsonTimeout, 403
∼CRoboDeviceNet, 39	94	RoboError_GilsonWrondID, 403
Axes_I, 401		RoboError_McsBus_UnknownCommand, 403
Axes_X, 401		RoboError_NoEndSwitch, 403
Axes_Y, 401		RoboError_NoMoreData, 403
Axes_Z, 401		RoboError_NoReference, 403
Axis_I, 401		RoboError_NoSpeedOrAcceleration, 403
Axis_X, 401		RoboError_OverPressure, 404
Axis_Y, 401		RoboError_ParameterNotAllowed, 404
Axis_Z, 401		RoboError PeristalticTimeout, 404
CancelPoolLoop, 394		RoboError_Phase0OutOfRange, 404
CancelPoolLoopAndS	topMovement, 394	RoboError_PollLoopCanceled, 404
CRoboDeviceNet, 394		RoboError_PollLoopCanceledAndStopMovement,
FindReference, 394		404
GetAirpressure, 395		RoboError_Pressure, 404
•	005	
GetAirpressureLimit, 3	995	RoboError_RangeExceeded, 404
GetAirValve, 395	05	RoboError_StateChangeNotPossible, 404
GetCurrentAirvalve, 39		RoboError_Timeout, 404
GetCurrentAirvalveLim		RoboError_UnknownCommand, 405
GetCurrentPosition, 39		RoboMainLowLevelCommand, 405
GetErrorAirpressure, 3		RoboStatusEvent, 405
GetErrorCurrentAirval		SetAirpressureLimit, 398
GetErrorVoltage12V, 3	396	SetAirValve, 399

SetCurrentAirvalveLimit, 399	IsPumpMotorOn, 407
SetCurrentAndAir, 399	m_pMcsBus_MotorControlNet, 408
SetInMovement, 399	m_pRoboFluidDevice, 408
SetMinPressure, 399	McsBus_MotorControl, 408
SetVoltage12VLimit, 399	PumpOff, 407
SetVoltage5VLimit, 399	PumpOn, 407
SetVoltageAirvalveLimit, 400	SetPumpSpeed, 407
SetVoltageRs485ALimit, 400	SetSingleValve, 408
SetVoltageRs485BLimit, 400	SetValve, 408
SetVoltageValvesLimit, 400	CRobolnjectDeviceNet, 409
StopMovement, 400	CRobolnjectDeviceNet, 409
CRoboDeviceNet::RoboMainLowLevelCommands, 608	CRoboocyte2DeviceNet, 409
FindReferencePhase0, 609	CRoboocyte2DeviceNet, 410
GetAxisConfig, 609	GetAxisLED, 410
GetHWConfig, 609	GetGilsonDevice, 410
GetHWRevision, 609	GetMcsBus Extension, 410
GetMaxNoPressure, 609	GetRoboDacg, 410
GetMaxNoPressureWaitTime, 609	GetRoboFluidDevice, 410
GetMaxPressureWaitTime, 609	SetAxisLED, 411
GetMinNoPressureWaitTime, 609	CRoboStatorDeviceNet, 411
GetMinPressure, 609	CRoboStatorDeviceNet, 412
GetMinPressureWaitTime, 610	FindReferencel, 412
GetParameter, 610	FindReferenceXY, 412, 413
GetPhases, 610	FindReferenceZ, 413
GetSearchReferenceFastAccel, 610	GetCurrentPositionI, 413
GetSearchReferenceFastSpeed, 610	GetCurrentPositionXY, 413
GetSearchReferenceFineAccel, 610	
,	GetCurrentPositionZ, 413
GetSearchReferenceFineSpeed, 610	HasRefl, 413
GetSearchReferenceMethod, 611	HasRefXY, 413
GetSearchReferenceMoveOut, 611	HasRefZ, 413
GetSearchReferenceOffsetPos, 611	MoveAbsI, 413, 414
GetUserParameter, 611	MoveAbsXY, 414
HasRef, 612	MoveAbsZ, 414
SetAxisConfig, 612	RoboMainStatorLowLevelCommand, 417
SetHWConfig, 612	SetAccelerationI, 414
SetHWRevision, 612	SetAccelerationNativeI, 414
SetMaxNoPressure, 612	SetAccelerationNativeXY, 414
SetMaxNoPressureWaitTime, 612	SetAccelerationNativeZ, 414
SetMaxPressureWaitTime, 612	SetAccelerationXY, 415
SetMinNoPressureWaitTime, 612	SetAccelerationZ, 415
SetMinPressure, 612	SetCurrentAndAirXY, 415
SetMinPressureWaitTime, 612	SetSpeedl, 415
SetParameter, 613	SetSpeedNativeI, 415
SetSearchReferenceFastAccel, 613	SetSpeedNativeXY, 415
SetSearchReferenceFastSpeed, 613	SetSpeedNativeZ, 415
SetSearchReferenceFineAccel, 613	SetSpeedXY, 415
SetSearchReferenceFineSpeed, 613	SetSpeedZ, 416
SetSearchReferenceMethod, 613	SetVelocityI, 416
SetSearchReferenceMoveOut, 614	SetVelocityXY, 416
SetSearchReferenceOffsetPos, 614	SetVelocityZ, 416
SetUserParameter, 614	StopMovementI, 416
CRoboFluidDeviceNet, 405	StopMovementXY, 416
\sim CRoboFluidDeviceNet, 406	StopMovementZ, 416
CloseAllValves, 406	CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands
CRoboFluidDeviceNet, 406	615
GetPumpSpeed, 407	FindReferencePhase0XY, 615
GetSingleValve, 407	CSafeISDeviceNet, 417
GetValve, 407	~CSafeISDeviceNet, 418

CSafeISDeviceNet, 418	SetAnalogOutChannels, 433
DacqDevice, 419	SetAnalogOutDACRange, 435
FluidControlDevice, 420	SetDacqLegacyMode, 435
RoboDevice, 420	SetReferenceElectrodeMode, 435
SetAdcChannels, 418	SetReferenceElectrodeSwitchState, 435
SetAdcSamplePos, 418	CSerialPortNet, 436
SetDacMode, 418	CSerialPortNet, 436
SetDacPeriode, 419	GetBytesAvailable, 437
SetDacPulseform, 419	Receive, 437
SetSwitches, 419	ReceiveString, 437
CSCUDacqGroupChannelSelectionNet, 420	Send, 437
CSCUDacqGroupChannelSelectionNet, 420	CStg200xBasicNet, 438
CSCUFunctionNet, 421	~CStg200xBasicNet, 442
!CSCUFunctionNet, 423	GetAnalogRanges, 442
~CSCUFunctionNet, 423	GetAnalogResolution, 442
AutomaticAnalogOut, 423	GetAutocalibrationDisabled, 443
CSCUFunctionNet, 423	GetAvailableMemory, 443
EnableAnalogOut, 424	GetBlankingEnable, 443, 444
GetAnalogOutADCRange, 424	GetCurrentRangeInNanoAmp, 444
GetAnalogOutChannels, 424	GetCurrentResolutionInNanoAmp, 444
GetAnalogOutDACRange, 424	GetDacAmplificationFactor, 445
GetAvailableHeadstages, 424	GetDACResolution, 445
GetAvailableHeadstagesEvent, 436	GetDiginValue, 445
GetFilterProperties, 425	GetDigoutMode, 445
GetFilterProperty, 425	GetDigoutValue, 446
GetHeadstageAdcBits, 425	GetElectrodeDacMux, 446
GetHeadstageAdcRangeInMicroVolt, 426	GetElectrodeEnable, 447
GetHeadstageDacBits, 426	GetElectrodeMode, 448
GetHeadstageDacCurrentRangeInMicroAmpere,	GetEnableAmplifierProtectionSwitch, 448, 449
426	GetExternalElectrodeEnable, 449
GetHeadstageDacCurrentResolutionInNanoAm-	GetFAAmplification, 450
pere, 427	GetHeadstage, 450
GetHeadstageDacVoltageRangeInMilliVolt, 427	GetListmodeIndexRange, 450
GetHeadstageDacVoltageResolutionInMicroVolt,	GetListmodeTriggerSource, 450
427	GetNumberOfAnalogChannels, 450
GetHeadstageGainInPermille, 428	GetNumberOfHWDACPaths, 450
GetHeadstageID, 428	GetNumberOfStimulationElectrodes, 450
GetHeadstageNumberOfAnalogChannels, 428	GetNumberOfStimulationSourcesPerElectrode,
GetHeadstageNumberOfStimulationChannels, 429	451
GetHeadstageSamplerate, 429	GetNumberOfSyncoutChannels, 451
GetHeadstageSerialNumber, 429	GetNumberOfTriggerInputs, 451
GetMaxNumberOfHeadstages, 430	GetOutputRate, 451
GetMaxStimulusChannelsPerHeadstage, 430	GetStgProgramInfo, 451, 452
GetReferenceElectrodeMode, 430	GetStgVersionInfo, 452
GetReferenceElectrodeSwitchState, 430	GetSyncoutMap, 452
HasAnalogOut, 431	GetTotalMemory, 453
HasGalvanicIsolation, 431	GetTriggerSource, 453
HasHSPowerSwitch, 431	GetVoltageRangeInMicroVolt, 453
IsAnalogOutEnabled, 431	GetVoltageResolutionInMicroVolt, 453
IsAutomaticAnalogOut, 432	ListModeSendStart, 454
IsHeadstageAvailable, 432	ListModeSendStop, 454
IsHeadstageAvailableEvent, 436	SendStart, 454
IsHSPowered, 432	SendStop, 454
IsInDacqLegacyMode, 432	SetAutocalibrationDisabled, 454
OnGetAvailableHeadstages, 433	SetBlankingEnable, 455, 456
OnIsHeadstageAvailable, 433	SetCurrentMode, 456
PowerHS, 433	SetDacAmplificationFactor, 456
SetAnalogOutADCRange, 433	SetDigoutMode, 457

	SetDigoutValue, 457		GetAvailableMemory, 487
	SetElectrodeDacMux, 457–459		GetCurrentRangeInNanoAmp, 488
	SetElectrodeEnable, 459–461		GetCurrentResolutionInNanoAmp, 488
	SetElectrodeMode, 462, 463		GetDACResolution, 488
	SetEnableAmplifierProtectionSwitch, 463, 464		GetMultiplexedDataChannelsInBlock, 488
	SetExternalElectrodeEnable, 465		GetNumberOfAnalogChannels, 489
	SetFAAmplification, 466		GetTotalMemory, 489
	SetHeadstage, 466		GetVoltageRangeInMicroVolt, 489
	SetListmodeIndexRange, 466		GetVoltageResolutionInMicroVolt, 489
	SetListmodeTriggerSource, 466		PollStatusEvent, 495
	SetMeasurementMode, 466		PrepareAndAppendData, 490
	SetOutputRate, 467		PrepareAndSendData, 491
	SetStgProgramInfo, 467		PrepareData, 491
	SetSyncoutMap, 467		SendMultiplexedData, 492
	SetTriggerSource, 467, 468		SendPreparedData, 492
	SetVoltageMode, 468		SendStart, 492
CSt	g200xDownloadBasicNet, 468		SendStop, 492
OOL	ClearChannelData, 470		SetupTrigger, 493
	ClearSyncData, 470		SetupTriggerSingle, 493
	DisableAutoReset, 470		StartPoll, 495
	EnableAutoReset, 471		StopPoll, 495
	ForceStatusEvent, 471	CSt	imulusFunctionNet::SidebandData, 617
	GetMemoryUsageDAC, 471	OSI	!SidebandData, 617
	GetMemoryUsageSyncout, 471		~SidebandData, 617
	GetSweepCount, 472		
	•		Duration, 617
	GetTrigger, 472		Sideband Peter 617
	ResetStatus, 472	CC+	SidebandData, 617
	SendChannelData, 474	USI	mulusFunctionNet::StimulusDeviceDataAndUnrolledData,
	SendSyncData, 474		618
	SetupRetriggerMode, 475		!StimulusDeviceDataAndUnrolledData, 619
	SetupTrigger, 475		~StimulusDeviceDataAndUnrolledData, 619
	SetupTriggerSingle, 476		Device Datal coath 610
004	Stimulus, 477		DeviceDataLength, 619
CSI	g200xDownloadNet, 477		StimulusDeviceDataAndUnrolledData, 619
	~CStg200xDownloadNet, 478		UnrolledAmplitude, 619
	ClearChannel_PrepareAndSendData, 479		UnrolledDuration, 619
	CStg200xDownloadNet, 478	00.	UnrolledSync, 619
	DisableMultiFileMode, 479	CSV	v2to64DeviceNet, 495
	EnableMultiFileMode, 479		~CSw2to64DeviceNet, 496
	GetModuleCurrent, 480		CSw2to64DeviceNet, 496
	GetModuleTemp, 480		GetChannel, 496
	MwPollStatusEvent, 483		GetChannels, 497
	PrepareAndAppendData, 480		GetNumber, 497
	PrepareAndSendData, 481		SetChannel, 497
	QueryTriggerstatus, 482	О.Т.	SetChannels, 497
	SendSegmentDefine, 482	Cic	xDeviceNet, 498
	SendSegmentSelect, 482		~CTcxDeviceNet, 500
	SendSegmentStart, 483		CalibrateThermocouple, 500
	SetOutputMap, 483		CTcxDeviceNet, 500
~~.	Stg200xPollStatusEvent, 483		FactoryReset, 500
CSt	imulusFunctionNet, 484		GetBoardTemp, 501
	ClearChannel_PrepareAndSendData, 486		GetCalibration, 501
	ClearChannelData, 486		GetCalibrationDecp, 501
	ClearMultiplexedData, 486		GetCalibrationMax, 501
	ClearSyncData, 486		GetCalibrationMin, 501
	CreateSideband, 486		GetCurrent, 501
	CStimulusFunctionNet, 485		GetD, 501
	ForceStatusEvent, 487		GetDDecp. 501

GetDevice, 501	SetEnableThermocouple, 509
GetDeviceType, 502	SetHeaterLimit, 509
GetDevname, 502	SetI, 509
GetDMax, 502	SetMaxHeaterPowerMultiwell, 510
GetDMin, 502	SetMaxP, 510
GetDuty, 502	SetOnOff, 510
GetEnableHeaterLimit, 502	SetP, 510
GetEnableThermocouple, 502	SetSensorType, 510
GetHasThermocouple, 502	SetSetpoint, 510
GetHeaterLimit, 503	SetThermocoupleNanovoltPerKelvin, 511
GetHeaterTemp, 503	CTEERFunctionNet, 511
Getl, 503	!CTEERFunctionNet, 513
GetIDecp, 503	~CTEERFunctionNet, 513
GetIMax, 503	CancelInternalCalibration, 513
GetlMin, 503	CTEERFunctionNet, 513
GetlOut, 503	GetAdapterCode, 514
GetMaxHeaterPowerMultiwell, 504	GetAdcOffsetU1, 514
GetMaxP, 504	GetAdcOffsetU2, 514
GetMaxpDecp, 504	GetAmplitude_nA, 514
GetMaxpBecp, 504 GetMaxpMax, 504	GetBytesPerSample, 514
GetMaxpMin, 504	GetClampMode, 514
GetNumControlChannels, 504	GetControllerParams, 515
GetNumDevices, 504	GetCurrentEnable, 515
GetNumMeasureChannels, 504	GetDacZero, 515
GetOnOff, 505	GetFrameErrorCounter, 515
GetP, 505	GetMayChunkSize Bute 516
GetPDecp, 505	GetMaxChunkSize_Byte, 516
GetPMax, 505	GetNumberOfAvailableSamples, 516
GetPMin, 505 GetPOut, 505	GetPeriod_us, 516 GetRotaryPositionCode, 516
GetPwrOut, 505	GetSampleBufferChunk, 516
GetPwrSet, 506	GetSampleBatter, 517
GetRes1, 506	GetSampleVoltageBuffer uV, 517
GetRes2, 506	GetScaleFactorU1, 517
GetResS, 506	GetScaleFactorU2, 517
GetResX, 506	GetUptimeSeconds, 518
GetROut, 506	GetWaveform, 518
GetSensorType, 506	IsInternalCalibrationFinished, 518
GetSetpoint, 506	IsSamplingFinished, 518
GetSetpointDecp, 507	SetAmplitude_nA, 518
GetSetpointMax, 507	SetBufferIndex, 519
GetSetpointMin, 507	SetClampMode, 519
GetThermocoupleCalibration, 507	SetControllerParams, 519
GetThermocoupleNanovoltPerKelvin, 507	SetCurrentEnable, 519
GetThermocoupleReferenceTemp, 507	SetExternalLED, 520
GetThermocoupleTemp, 507	SetLiquidResistance, 520
GetThermocoupleTempAbs, 508	SetPeriod us, 520
GetUnit, 508	SetWaveform, 520
GetUOut, 508	StartInternalCalibration, 520
GetValue, 508	StartSampling, 521
GetValueHires, 508	StopSampling, 521
GetVolti, 508	CTEERMachineDeviceNet, 521
SetCalibration, 508	~CTEERMachineDeviceNet, 522
SetD, 509	CTEERMachineDeviceNet, 521
SetDevice, 509	TEERFunctionNet, 522
SetDeviceType, 509	CurrentRangeInNanoAmp
SetDevname, 509	W2100_StimulusParametersNet, 621
SetEnableHeaterLimit, 509	CurrentResolutionInNanoAmp
	•

W2100_StimulusParametersNet, 621	amplification, 29
CUsbExceptionNet, 522	channel, 29
CUsbExceptionNet, 522, 523	dacqgroup, 29
Status, 523	CW2100_StimulatorFunctionNet, 531
CutoffFrequency	BOOST_BIT, 538
CCreateFilterNet, 49	ClearChannelData, 533
CVoltageRangeInfoNet	CW2100_StimulatorFunctionNet, 533
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	
523	GetBoostPreTime, 533
CW2100_FunctionNet, 524	GetCurrentRangeInNanoAmp, 533
ClearStimulusParametersCache, 526	GetCurrentResolutionInNanoAmp, 534
ClearUserDefinedNameCache, 526	GetDACResolution, 534
CW2100_FunctionNet, 525	GetDigitalStimulatorTrigger, 534
DeselectAllHeadstages, 526	GetDigitalStimulatorTriggerSlope, 534
DeselectHeadstage, 526	GetNumberOfAnalogChannels, 534
GetAccelGyroCurrentRate, 526	GetNumberOfSyncoutChannels, 535
GetAccelGyroDesiredRate, 526	GetNumberOfTriggerInputs, 535
GetAccelGyroEnabled, 526	GetStimulationPatternMemory, 535
GetAccelRange, 526	GetTimeResolutionInNanoSeconds, 535
GetAnalogOutChannel, 527	GetTimeSlot, 535
GetAnalogOutFilter, 527	GetVoltageRangeInMicroVolt, 535
GetAudioChannels, 527	GetVoltageResolutionInMicroVolt, 535
GetAvailableHeadstages, 527	GND_SWITCH_BIT, 539
GetBatteryState, 527	PollStatusEvent, 539
GetDacRange, 527	PrepareData, 537
GetFilterProperties, 527	PrepareDataSync, 537
GetFilterProperty, 527	SelectTimeSlot, 537
GetFPGAFirmwareType, 527	SendPreparedData, 537
GetGyroRange, 527	SendStart, 537
GetHeadstageOnOff, 528	SendStop, 538
GetHeadstageSamplingActive, 528	SetBoostAlwaysOnMode, 538
GetMultiHeadstageMode, 528	SetDigitalStimulatorTrigger, 538
GetPicFirmwareType, 528	SetDigitalStimulatorTriggerSlope, 538
GetSelectedChannels, 528	StartPoll, 538
GetSelectedHeadstageState, 528	StopPoll, 538
GetStimulusParametersCache, 528	SYNC_BIT0, 539
GetStimulusParametersFromSelectedHS, 528	SYNC_BIT1, 539
GetStiumlusParameters, 529	CW2100DacqGroupChannelSelectionNet, 539
GetUserDefinedName, 529	CW2100DacqGroupChannelSelectionNet, 539
GetUserDefinedNameCache, 529	CWarnerUssingDeviceNet, 540
GetUserDefinedNameFromSelectedHS, 529	!CWarnerUssingDeviceNet, 540
PulseGenerator, 531	~CWarnerUssingDeviceNet, 540
SelectHeadstage, 529	CWarnerUssingDeviceNet, 540
SetAccelGyroDesiredRate, 529	WarnerUssingFunction, 541
SetAccelGyroEnabled, 529	CWarnerUssingFunctionNet, 541
SetAccelRange, 529	!CWarnerUssingFunctionNet, 543
SetAnalogOutChannel, 530	~CWarnerUssingFunctionNet, 543
SetAnalogOutFilter, 530	CompensateElectrodeOffset, 543
SetAudioChannels, 530	CWarnerUssingFunctionNet, 543
SetDacRange, 530	GetAvailableChambers, 544
SetGyroRange, 530	GetChannelsCountOfChamber, 544
SetHeadstageOnOff, 530	GetClampMode, 544
SetHeadstageSamplingActive, 530	GetComplianceVoltage, 544
SetHeadstageToSleep, 530	GetDacNampsPerDigitHighCurrentRange, 545
SetMultiHeadstageMode, 531	GetDacPampsPerDigitLowCurrentRange, 545
SetSelectedChannels, 531	GetDacZero, 545
Stimulator, 531	GetHighCurrentRange, 546
CW2100_FunctionNet::AudioChannelsNet, 29	GetIdleModeOffset, 546

GetLiquidResistance, 546	GetTotalNumberOfDigitalPorts, 564
GetLowCurrentRange, 547	GetTotalNumberOfTables, 564
GetNumberOfAvailableChambers, 547	GetTotalNumberOfValves, 564
GetNumberOfHardwareSlotsForChambers, 547	GetTotalTableSize, 564
GetU1Offset, 547	GetValveActive, 564
GetU1Reference, 548	GetValveActiveEvent, 578
GetU2Offset, 548	GetValveBoardRevision, 565
GetU2Reference, 548	GetValveBoardRevisionEvent, 578
GetUnitDescription, 549	GetValveBoardRevisionString, 565
GetUnitExponent, 549	GetValveDigitalInPort, 565
GetUnitName, 549	GetValveDigitaIInPortEvent, 578
GetUnitsPerDigit, 550	GetValveLedOn, 565
GetUptimeSeconds, 550	GetValveLedOnEvent, 579
GetVoltageClampControllerParam_D, 550	GetValveManualGroup, 566
GetVoltageClampControllerParam_I, 551	GetValveManualGroupEvent, 579
GetVoltageClampControllerParam_P, 551	GetValveManualState, 566
IsChamberAvailable, 551	GetValveManualStateEvent, 579
IsHighCurrentMode, 552	GetValveMode, 566
IsInternalCalibrationFinished, 552	GetValveModeEvent, 579
IsPulseEnabled, 552	GetValvesActiveMap, 566
SetClampMode, 553	GetValvesManualStateMap, 567
SetEnablePulse, 553	GetValveTableEntry, 567
SetHighCurrentMode, 553	IsDigitalOutPortInverted, 567
SetIdleModeOffset, 554	IsDigitalOutPortInvertedEvent, 579
SetLiquidResistance, 554	IsValveDigitalInInverted, 567
SetLowCurrentMode, 554	IsValveDigitalInInvertedEvent, 579
SetPulse, 554	IsValveOpen, 568
SetVoltageClampControllerParam_D, 555	IsValveOpenEvent, 579
SetVoltageClampControllerParam_I, 555	IsValveOpenInAnalogMode, 568
SetVoltageClampControllerParam_P, 555	IsValveOpenInAnalogModeEvent, 579
WaitForAllChambers, 556	IsValveOpenInDigitalMode, 568
WaitForChamber, 556	IsValveOpenInDigitalModeEvent, 580
CWarnerValveControllerDeviceNet, 556	LoadValveTable, 569
!CWarnerValveControllerDeviceNet, 560	OnGetActiveRunningTableNumber, 569
~CWarnerValveControllerDeviceNet, 560	OnGetAnalogThresholdHigh, 569
ClearTableName, 561	OnGetAnalogThresholdLow, 569
ClearValveTable, 561	OnGetAnalogVoltage, 569
CWarnerValveControllerDeviceNet, 560	OnGetCurrentNumberOfValves, 569
GetActiveRunningTableNumber, 561	OnGetDigitalOutPortValve, 569
GetActiveRunningTableNumberEvent, 577	OnGetDigitalPortDirection, 569
GetAnalogThresholdHigh, 561	OnGetDisplayMode, 570
GetAnalogThresholdHighEvent, 577	OnGetTableNamebyIndex, 570
GetAnalogThresholdLow, 561	OnGetValveActive, 570
GetAnalogThresholdLowEvent, 577	OnGetValveBoardRevision, 570
GetAnalogVoltage, 562	OnGetValveDigitalInPort, 570
GetAnalogVoltageEvent, 577	OnGetValveLedOn, 570
GetCurrentEditTableNumber, 562	OnGetValveLedon, 570 OnGetValveManualGroup, 570
GetCurrentNumberOfValves, 562	OnGetValveManualState, 570
GetCurrentNumberOfValves, 362 GetCurrentNumberOfValvesEvent, 577	OnGetValveMode, 570
GetDigitalOutPortValve, 562	OnlsDigitalOutPortInverted, 571
_	-
GetDigitalOutPortValveEvent, 578	OnlsValveOpen 571
GetDigitalPortDirection, 563	OnlsValveOpenInAppledMode, 571
GetDigitalPortDirectionEvent, 578	OnlsValveOpenInAnalogMode, 571
GetDisplayMode, 563	OnlsValveOpenInDigitalMode, 571
GetDisplayModeEvent, 578	OnTableEntryChanged, 571
GetTableName, 563	SetActiveRunningTableNumber, 571
GetTableNamebyIndex, 563	SetAnalogThresholdHigh, 572
GetTableNamebyIndexEvent, 578	SetAnalogThresholdLow, 572

SetCurrentEditTableNumber, 572	SetResetFilter, 586
SetDefault, 572	SetRFFrequencyHeadstage, 586
SetDigitalOutPortInvert, 572	SetRFFrequencyReceiver, 587
SetDigitalOutPortValve, 573	SetRFFrequencyReceiverEeprom, 587
SetDigitalPortDirection, 573	SetRFLostBehaviour, 587
SetDisplayMode, 573	SetRFPower, 587
SetTableName, 573	SetSelectedHeadstage, 587
SetTableStep, 574	SetSerialNumberHeadstage, 587
SetTableStepAll, 574	SetWPADebugMode, 587
SetValveActive, 574	SetWPAType, 587
SetValveDigitalInInvert, 574	CWirelessBaseFunctionNet, 588
SetValveDigitalInPort, 575	CreateWirelessHeadstageSerialNumberString,
SetValveLedOn, 575	588
SetValveManualGroup, 575	CWirelessBaseFunctionNet, 588
SetValveManualState, 575	Ovviiciossbasei unotionivet, soo
SetValveMode, 576	DacqDevice
SetValvesActiveMap, 576	CSafeISDeviceNet, 419
SetValvesManualStateMap, 576	dacqgroup
·	CW2100_FunctionNet::AudioChannelsNet, 29
SetValveTableEntry, 576	DACResolution
StoreValveTable, 577	W2100_StimulusParametersNet, 621
TableEntryChangedEvent, 580 CWarnerValveControllerDeviceTesterFunctionNet, 580	DataState
•	HeadStageIDTypeState, 604
!CWarnerValveControllerDeviceTesterFunctionNet,	DeepCopy
581	CCMCCMap Davida Nature Davida Officta ve at David
~CWarnerValveControllerDeviceTesterFunctionNet,	366
581	DefineAmplification
CWarnerValveControllerDeviceTesterFunctionNet,	CPgaDeviceNet, 331
581	DefineFrequencyRange
GetIO, 581	CPgaDeviceNet, 331
GetSync, 581	DefineNumAmplifications
SetADC, 582	CPgaDeviceNet, 331
SetIO, 582	DefineNumFrequencyRanges
SetIODirection, 582	CPgaDeviceNet, 331
SetTrigger, 582	DeselectAllHeadstages
SetTriggerSyncDirection, 583	CW2100_FunctionNet, 526
CWClassicFunctionNet, 583	
CWClassicFunctionNet, 584	DeselectHeadstage CW2100 FunctionNet, 526
GetFilterParametersHeadstage, 584	
GetHasChecksum, 584	DetectChipType
GetHasRedLedHeadstage, 584	CCMOSMea_FunctionNet, 36
GetHeadstageOnOff, 584	DEVICE_NOT_FOUND
GetResetFilter, 585	Mcs::Usb, 27
GetRFConnectionStatus, 585	DeviceArrival
GetRFFrequencyHeadstage, 585	CMcsUsbListNet, 239
GetRFFrequencyReceiver, 585	DeviceData
GetRFPower, 585	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData
GetScanHeadstagesResult, 585	619
GetSelectedHeadstage, 585	DeviceDataLength
GetSerialNumberHeadstage, 585	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData
GetWPADebugMode, 585	619
GetWPAType, 585	DeviceId
ResetChannelmap, 586	CMcsUsbListEntryNet, 236
ScanForHeadstages, 586	DeviceIdNet, 588
SetChannelmap, 586	BcdDevice, 590
SetFilterParametersHeadstage, 586	BusType, 590
SetHasChecksum, 586	DeviceIdNet, 589
SetHeadstageOnOff, 586	IdProduct, 590
SetHWSelectedChannels, 586	IdVendor, 590
Gett tyvoeteotedottatillets, Jou	operator- 580

DeviceName	CRoboDacqNet, 381
CMcsUsbListEntryNet, 236	Emu_SetCellPotential
DeviceRemoval	CRoboDacqNet, 381
CMcsUsbListNet, 239	Emu_SetCellResists
DigitalSource	CRoboDacqNet, 381
DigitalSource< digitalsourceenum >, 590	Emu_SetElectrodeResists
DigitalSource< digitalsourceenum >, 590	CRoboDacqNet, 381
DigitalSource, 590	Emu_SetNoise
MaxBitNumber, 590	CRoboDacqNet, 381
size, 591	EnableAnalogOut
Source, 591	CSCUFunctionNet, 424
DisableAutoReset	EnableAutoReset
CStg200xDownloadBasicNet, 470	CStg200xDownloadBasicNet, 471
DisableMultiFileMode	EnableChannelsInGroup
CStg200xDownloadNet, 479	CCMOSMea_FunctionNet, 36, 37
Disconnect	CDacqGroupChannelSelectionTemplateNet< Dac
CMcsUsbNet, 246	qGroupChannelEnumTemplateNet, Dac
DisConnectDevice	qGroupChannelEnumTemplate, CDevice
CRadioControledDevicesNet, 365	GroupChannelInfoTemplateNet >, 54
DoRamp	EnableChecksum
CRoboDacqNet, 380	CMeaDeviceNet, 278
DownloadFirmware	COctoPotDeviceNet, 319
CMcsUsbFactoryNet, 224	EnableDigitaIIn
DriverVersionNet, 591	CMeaDeviceNet, 278, 279
∼DriverVersionNet, 592	COctoPotDeviceNet, 319
DriverVersionNet, 592	EnableExceptions
DriverVersionNet::FormatVersion, 592	CMcsUsbNet, 246
GetDestinationCode, 592	EnableMultiFileMode
GetDestinationName, 592, 594	CStg200xDownloadNet, 479
GetMajor, 594	EnableTimestamp
GetMinor, 594, 595	CMeaDeviceNet, 279
GetNumEntries, 595	COctoPotDeviceNet, 319
GetStatus, 595	EnableUserTrigger
GetVersionInt, 595, 596	CLIH3DeviceNet, 122
GetVersionString, 596	enCMosMeaChipType
DriverVersionNet::FormatVersion	Mcs::Usb, 27
DriverVersionNet, 592	EnSTG200x_STATUS
DSP	Mcs::Usb, 27
FirmwareDestinationNames, 597	Entry
DummyCommand	HeadStageIDType, 601
CLIH3DeviceNet, 122	Equals
Duration	CMcsUsbListEntryNet, 233
CStimulusFunctionNet::SidebandData, 617	HeadStageIDType, 601
,	HeadstageIDTypeObject, 603
ElectricalStimulation	EraseEepromRegisterPreconfig
HeadStageIDType, 601	CMcsUsbNet, 246
EmptyKey	EraseFilterParameterPermanent
CMcsUsbNet, 246	CFilterConfigurationNet, 69
Emu_GetCellCapacity	CFilterConfigurationRegisterNet, 70, 71
CRoboDacqNet, 380	ErasePermanentAdcOffset
Emu_GetCellPotential	CLIH3DeviceNet, 123
CRoboDacqNet, 380	ErasePermanentDacOffset
Emu_GetCellResists	CLIH3DeviceNet, 123
CRoboDacqNet, 380	Error_Callback_Aquisition_Stopped
Emu_GetElectrodeResists	CMcsUsbDacqNet, 218
CRoboDacqNet, 381	Error_Callback_Data_lost
Emu_GetNoise	CMcsUsbDacqNet, 218
CRoboDacqNet, 381	Error_Callback_Frames_Lost
Emu_SetCellCapacity	

CMcsUsbDacqNet, 218	CFilterPropertyNet, 73
Error_Callback_Packet_Error	FilterType
CMcsUsbDacqNet, 218	CFilterPropertyNet, 73
Error_Callback_Queue_Full	FindFilter
CMcsUsbDacqNet, 219	CCreateFilterNet, 49
Error_Callback_RingQueue_Full	FindFirmwareVersionMagicInBuffer
CMcsUsbDacqNet, 219	CMcsUsbFactoryNet, 224
ErrorEvent	FindReference
CMcsUsbDacqNet, 219	CRoboDeviceNet, 394
•	FindReferencel
FactoryReset	CRoboStatorDeviceNet, 412
CTcxDeviceNet, 500	FindReferencePhase0
FeedbackGetSampleTimerCount	CRoboDeviceNet::RoboMainLowLevelCommands,
CMeaFeedbackFunctionNet, 287	609
FeedbackSetAnalogSource	FindReferencePhase0XY
CMeaFeedbackFunctionNet, 287	CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands,
FeedbackSetChannelFilter	615
CMeaFeedbackFunctionNet, 287	FindReferenceXY
FeedbackSetDigitalMapping	CRoboStatorDeviceNet, 412, 413
CMeaFeedbackFunctionNet, 287	FindReferenceZ
FeedbackSetFeedback	CRoboStatorDeviceNet, 413
CMeaFeedbackFunctionNet, 287	FirePressurePulse
FeedbackSetFilterOff	CPPCFunctionNet, 347
CMeaFeedbackFunctionNet, 288	FirmwareDestinationNames, 596
FeedbackSetFilterParameter	Altera, 597
CMeaFeedbackFunctionNet, 288	Bootstrap, 597
FeedbackSetFilterParameter32	BUS1_MCSBUS1, 597
CMeaFeedbackFunctionNet, 288	BUS1_MCSBUS2, 597
FeedbackSetGlobalChannelFilter	DSP, 597
CMeaFeedbackFunctionNet, 288	FPGA2, 597
FeedbackSetIIRFilterParameter	FPGA3, 598
CMeaFeedbackFunctionNet, 288	FPGA4, 598
FeedbackSetLogic	FPGA5, 598
CMeaFeedbackFunctionNet, 288	FPGA6, 598
FeedbackSetMkFilter	MCSBUS1, 598
CMeaFeedbackFunctionNet, 288	MCSBUS10, 598
FeedbackSetNumberOfLogics	MCSBUS11, 598
CMeaFeedbackFunctionNet, 289	MCSBUS12, 598
FeedbackSetNumberOfRateCounter	MCSBUS13, 598
CMeaFeedbackFunctionNet, 289	MCSBUS2, 598
FeedbackSetNumberOfRateDetectors	MCSBUS3, 598
CMeaFeedbackFunctionNet, 289	MCSBUS4, 599
FeedbackSetNumberOfSpikeDetectors	MCSBUS5, 599
CMeaFeedbackFunctionNet, 289	MCSBUS6, 599
FeedbackSetNumberOfTriggers	MCSBUS7, 599
CMeaFeedbackFunctionNet, 289	MCSBUS8, 599
FeedbackSetRateCounter	MCSBUS9, 599
CMeaFeedbackFunctionNet, 289	MCU1, 599
FeedbackSetRateDetector	PIC, 599
CMeaFeedbackFunctionNet, 289	PIC2, 599
FeedbackSetSpikeDetectorThreshold	PIC3, 599
CMeaFeedbackFunctionNet, 289	PIC4, 599
FeedbackSetTrigger	USB, 600
CMeaFeedbackFunctionNet, 290	FluidControlDevice
FilterActive	CSafeISDeviceNet, 420
CFilterPropertyNet, 73	ForceStatusEvent
FilterBand	CStg200xDownloadBasicNet, 471
CFilterPropertyNet, 73	CStimulusFunctionNet, 487
FilterFamily	

FPGA2	GetAdapterCode
FirmwareDestinationNames, 597	CMealmpedanceDeviceNet, 291
FPGA3	CTEERFunctionNet, 514
FirmwareDestinationNames, 598	GetAdapterType
FPGA4	CMcsUsbDacqNet, 194
FirmwareDestinationNames, 598	GetAdc
FPGA5	CFluidControlDeviceNet, 75
	GetAdcDataFormat
FirmwareDestinationNames, 598	
FPGA6	CMcsUsbDacqNet, 194
FirmwareDestinationNames, 598	GetADCInputOffset
FromIntPtr	CCMOSMea_FunctionNet, 37
StgStatusNet, 618	GetAdcOffset
FromPtr	CLIH3DeviceNet, 123
StgStatusNet, 618	COctoPotDeviceNet, 319
FX3MCSDataAddress	GetAdcOffsetU1
CMcsUsbFactoryNet, 229	CTEERFunctionNet, 514
FX3MCSDataDeviceIdOffset	GetAdcOffsetU2
CMcsUsbFactoryNet, 229	CTEERFunctionNet, 514
FX3MCSDataIFB1ImageOffset	GetADCs
CMcsUsbFactoryNet, 229	CMcsBus SensorNet, 155
FX3MCSDataIFB2ImageOffset	GetADCsLoop
CMcsUsbFactoryNet, 229	CMcsBus_SensorNet, 155
FX3MCSDataVersionOffset	GetAdcZero
CMcsUsbFactoryNet, 229	CMcsUsbDacqNet, 194
FYIProgram	GetAirpressure
	CRoboDeviceNet, 395
CFYIDeviceNet, 81	
FYITemp	GetAirpressureLimit
CFYIDeviceNet, 81	CRoboDeviceNet, 395
Gain	GetAirValve
CMeaDeviceNet, 283	CRoboDeviceNet, 395
	GetAllDigout
Get2AnalogInput	CRoboDacqNet, 381
CMcsBus_SensorNet, 154	GetAmplification
Get2DigitalInput	CPgaDeviceNet, 331
CMcsBus_SensorNet, 154	GetAmplitude_nA
Get4ADC	CTEERFunctionNet, 514
CMcsBus_SensorNet, 154	GetAnalogGain
Get4ADCAverage	CMeaDeviceNet, 280
CMcsBus_SensorNet, 154	GetAnalogOutADCRange
Get4ADCCatchampAverageShift	CSCUFunctionNet, 424
CMcsBus_SensorNet, 154	GetAnalogOutChannel
Get4ADCMode	CW2100 FunctionNet, 527
CMcsBus_SensorNet, 155	GetAnalogOutChannels
Get4DAC	CSCUFunctionNet, 424
CMcsBus_SensorNet, 155	GetAnalogOutDACRange
GetAbsMaxCurrentInMicroAmp	
CMultiwellOptoStimFunctionNet, 313	CSCUFunctionNet, 424
GetAccelGyroCurrentRate	GetAnalogOutFilter
CW2100 FunctionNet, 526	CW2100_FunctionNet, 527
GetAccelGyroDesiredRate	GetAnalogRanges
-	CStg200xBasicNet, 442
CW2100_FunctionNet, 526	GetAnalogResolution
GetAccelGyroEnabled	CStg200xBasicNet, 442
CW2100_FunctionNet, 526	GetAnalogThresholdHigh
GetAccelRange	CWarnerValveControllerDeviceNet, 561
CW2100_FunctionNet, 526	GetAnalogThresholdHighEvent
GetActiveRunningTableNumber	CWarnerValveControllerDeviceNet, 577
CWarnerValveControllerDeviceNet, 561	GetAnalogThresholdLow
GetActiveRunningTableNumberEvent	CWarnerValveControllerDeviceNet, 561
CWarnerValveControllerDeviceNet, 577	2 2 2 2 2 2 2 2 2 2 2 2 2 2 2 2 2 2 2 2

GetAnalogThresholdLowEvent	GetBaseFrequency
CWarnerValveControllerDeviceNet, 577	CRFFunctionNet, 372
GetAnalogValueUnit	GetBaseSamplerate
CMcsUsbDacqNet, 195	CCMOSMeaDeviceNet, 46
GetAnalogVoltage	GetBath
CPPCFunctionNet, 347	CCMOSMea_FunctionNet, 37
CPPS_FunctionNet, 356	GetBathMode
CWarnerValveControllerDeviceNet, 562	CCMOSMea_FunctionNet, 37
GetAnalogVoltageEvent	GetBatteryState
CWarnerValveControllerDeviceNet, 577	CW2100 FunctionNet, 527
GetAnalogVoltageRange	GetBatteryVoltage
CPPCFunctionNet, 347	CMultiBatteryChargerDeviceNet, 299
GetAnalogVoltages	GetBiQuad
CPPS_FunctionNet, 356	CCreateFilterNet, 49
GetArraySize	GetBiQuads
CMealmpedanceDeviceNet, 291	CCreateFilterNet, 49
GetAudioChannels	GetBlankingEnable
CMeaAudioFunctionNet, 264, 265	CStg200xBasicNet, 443, 444
CW2100 FunctionNet, 527	GetBoardTemp
GetAudioOutDacParameter	CTcxDeviceNet, 501
CLIH3DeviceNet, 123	GetBoostAlwaysOnMode
GetAutocalibrationDisabled	CW2100 StimulatorFunctionNet, 533
CStg200xBasicNet, 443	GetBoostPreTime
GetAvailableBaseSamplerates	CW2100_StimulatorFunctionNet, 533
CCMOSMeaDeviceNet, 46	GetBubbleState
GetAvailableChambers	
	CPPS_FunctionNet, 356
CWarnerUssingFunctionNet, 544	GetBubbleStatus
GetAvailableDeviceList	CMcsBus_SensorNet, 155
CRFFunctionNet, 371	GetBuffer
GetAvailableDeviceListEx	CGenericDevelopDeviceNet, 89
CRFFunctionNet, 371	GetBusAddress
GetAvailableHeadstages	CMcsBusNet, 169
CSCUFunctionNet, 424	GetBusAddressEeprom
CW2100_FunctionNet, 527	CMcsBusNet, 169
GetAvailableHeadstagesEvent	GetByteBuffer
CSCUFunctionNet, 436	CGenericDevelopDeviceNet, 89
GetAvailableMemory	GetBytesAvailable
CStg200xBasicNet, 443	CSerialPortNet, 437
CStimulusFunctionNet, 487	GetBytesPerSample
GetAvailableSampleRates	CTEERFunctionNet, 514
CMcsUsbDacqNet::CHWInfo, 113	GetCalibration
GetAvailableStateList	CTcxDeviceNet, 501
CRFFunctionNet, 371	GetCalibrationDecp
GetAvailableStateListEx	CTcxDeviceNet, 501
CRFFunctionNet, 371	GetCalibrationMax
GetAvailableVoltageRangesInMicroVolt	CTcxDeviceNet, 501
CMcsUsbDacqNet::CHWInfo, 113	GetCalibrationMin
GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVo	
CMcsUsbDacqNet::CHWInfo, 113	GetCapacityC
GetAxisConfig	CRoboDacqNet, 381
CRoboDeviceNet::RoboMainLowLevelCommands,	GetCapacityV
609	CRoboDacqNet, 381
GetAxisLED	GetCapacityX
CRoboocyte2DeviceNet, 410	CRoboDacqNet, 382
	·
GetAxisParametersSignedEeprom	GetCardinalDacqSamplerate
CMcsBus_AxisParametersNet, 130	ClnterfaceboardFunctionNet, 119
GetAxisParametersUnsignedEeprom	GetCardinalStgOutputrate
CMcsBus AxisParametersNet. 130	CInterfaceboardFunctionNet, 119

CotChannal	CotConfiguration Dit Ava
GetChannel	GetConfigurationBitAxc
CSw2to64DeviceNet, 496	CRoboDacqNet, 382
GetChannelDataFillSize	GetConfigurationBitBlu_Led
CMcsUsbDacqNet, 195	CRoboDacqNet, 382
GetChannelDatal16	GetConfigurationBitBlu_LedToggleFast
CCMOSMeaDeviceNet, 46	CRoboDacqNet, 382
GetChannelDatal32	GetConfigurationBitBlu_LedToggleSlow
CCMOSMeaDeviceNet, 46	CRoboDacqNet, 382
GetChannelDataUI16	GetConfigurationBitCC_Gen
CCMOSMeaDeviceNet, 46	CRoboDacqNet, 382
GetChannelDataUI32	GetConfigurationBitCV_Gen
CCMOSMeaDeviceNet, 47	CRoboDacqNet, 382
GetChannelLayout	GetConfigurationBitRC_Gen
CMcsUsbDacqNet, 195	CRoboDacqNet, 382
GetChannels	GetConfigurationBitRed_Led
CMultiBatteryChargerDeviceNet, 299	CRoboDacqNet, 383
CSw2to64DeviceNet, 497	GetConfigurationBitRed_LedSaturation
GetChannelsCountOfChamber	CRoboDacqNet, 383
CWarnerUssingFunctionNet, 544	GetConfigurationBitRed_LedToggleFast
GetChannelsInBlock	CRoboDacqNet, 383
CMcsUsbDacqNet, 195	GetConfigurationBitRed_LedToggleSlow
GetChannelState	CRoboDacqNet, 383
CMultiBatteryChargerDeviceNet, 299	GetConfigurationBitRelais
GetChargeCapacity	CRoboDacqNet, 383
CMultiBatteryChargerDeviceNet, 300	GetConfigurationBitRV_Gen
GetChargeCurrent	CRoboDacqNet, 383
CMultiBatteryChargerDeviceNet, 300	GetConfigurationBits
GetChargingMode	CRoboDacqNet, 383
CMultiBatteryChargerDeviceNet, 300	GetConfigurationBitStream
GetChargingPCoefficient	CRoboDacqNet, 383
CMultiBatteryChargerDeviceNet, 300	GetConfigurationBitSupply
GetChecksumFromFX3Image	CRoboDacqNet, 383
CMcsUsbFactoryNet, 224	GetConnectedDevice
GetCheckVoltage	CRFFunctionNet, 372
COkuvisionStimulatorDeviceNet, 323	GetControllerParams
GetClampAmpSerialNumber	CTEERFunctionNet, 515
CRoboDacqNet, 382	GetCrossTalkOffset
GetClampMode	CRoboDacqNet, 383
CTEERFunctionNet, 514	GetCrossTalkOptimum
CWarnerUssingFunctionNet, 544	CRoboDacqNet, 383
GetCMOSDataDictionary	GetCurrent
CCMOSMeaDeviceNet, 47	CTcxDeviceNet, 501
GetCoilCommunication	GetCurrentAirvalve
CPositionIIDeviceNet, 334	CRoboDeviceNet, 395
GetColorRgb	GetCurrentAirvalveLimit
CMultiwellOptoStimFunctionNet, 313	CRoboDeviceNet, 395
GetColorStr	GetCurrentCycle
CMultiwellOptoStimFunctionNet, 314	CMeaCoatDeviceNet, 272
GetCommand	GetCurrentEditTableNumber
CMcsBusNet, 169, 170	CWarnerValveControllerDeviceNet, 562
CPedoterDeviceNet, 328	GetCurrentEnable
CRoboDacqNet, 382	CTEERFunctionNet, 515
GetComplianceVoltage	GetCurrentFactor
CWarnerUssingFunctionNet, 544	COkuvisionStimulatorDeviceNet, 323
GetConfiguration	GetCurrentNumberOfValves
CMcsUsbNet, 247	CWarnerValveControllerDeviceNet, 562
GetConfigurationBit	GetCurrentNumberOfValvesEvent
CRoboDacqNet, 382	CWarnerValveControllerDeviceNet, 577

GetCurrentPosition	CTcxDeviceNet, 501
CRoboDeviceNet, 395	GetDebugData
GetCurrentPositionI	CPositionIIDeviceNet, 334
CRoboStatorDeviceNet, 413	GetDestination
GetCurrentPositionXY	CMcsUsbFactoryNet, 224
CRoboStatorDeviceNet, 413	GetDestinationCode
GetCurrentPositionZ	DriverVersionNet, 592
CRoboStatorDeviceNet, 413	GetDestinationDisplayLabel
GetCurrentRangeInNanoAmp	CMcsUsbFactoryNet, 224
CStg200xBasicNet, 444	GetDestinationName
CStimulusFunctionNet, 488	CMcsUsbFactoryNet, 224, 225
CW2100_StimulatorFunctionNet, 533	DriverVersionNet, 592, 594
GetCurrentResolutionInNanoAmp	GetDestinationSerialNumber
CStg200xBasicNet, 444	CMcsUsbFactoryNet, 225
CStimulusFunctionNet, 488	GetDestinationTargetAddress
CW2100_StimulatorFunctionNet, 534	CMcsUsbFactoryNet, 225
GetCycle	GetDetectionThreshold
CMeaCleanDeviceNet, 268	CMcsBus_SensorNet, 155
GetCycles	GetDetectorValue
CMeaCleanDeviceNet, 268	CMcsBus_SensorNet, 155
CMeaCoatDeviceNet, 272	GetDevice
GetD	CTcxDeviceNet, 501
CTcxDeviceNet, 501	GetDeviceCannotStallOutRequests
GetDacAmplificationFactor	CMcsUsbNet, 247
CStg200xBasicNet, 445	GetDeviceCapableSpeed
GetDacIdleValue	CMcsUsbNet, 247
CLIH3DeviceNet, 124	GetDeviceEnum
GetDacNampsPerDigitHighCurrentRange	CMcsUsbNet, 247
CWarnerUssingFunctionNet, 545	GetDeviceGroupChannelInfos
GetDACOffset	CDacqGroupChannelSelectionTemplateNet< Dac-
CGrapheneFunctionNet, 103, 104	qGroupChannelEnumTemplateNet, Dac-
COkuvisionStimulatorDeviceNet, 323	qGroupChannelEnumTemplate, CDevice-
GetDacOffset	GroupChannelInfoTemplateNet >, 54
CDacCalibrationFunctionNet, 51	GetDeviceId
CLIH3DeviceNet, 124	CMcsUsbNet, 247
COctoPotDeviceNet, 319	GetDeviceList
GetDacPampsPerDigitLowCurrentRange	CPositionImpDeviceNet, 342
CWarnerUssingFunctionNet, 545	GetDeviceNames
GetDacqRunStatus	CRadioControledDevicesNet, 365
CLIH3DeviceNet, 124	GetDeviceRootHubVendorEnum
GetDacRange	CMcsUsbNet, 247
CW2100 FunctionNet, 527	GetDeviceRootHubVendorID
GetDACResolution	CMcsUsbNet, 247
CStg200xBasicNet, 445	GetDeviceRootHubVendorName
CStimulusFunctionNet, 488	CMcsUsbNet, 247
CW2100 StimulatorFunctionNet, 534	GetDeviceSpeed
GetDACs	CMcsUsbNet, 248
CMcsBus_SensorNet, 155	GetDeviceType
GetDacUseIdleValue	CTcxDeviceNet, 502
CLIH3DeviceNet, 124	GetDevname
GetDacZero	CTcxDeviceNet, 502
CTEERFunctionNet, 515	GetDigin
CWarnerUssingFunctionNet, 545	CFluidControlDeviceNet, 75
GetDataFormat	GetDigInState
CMcsUsbDacqNet, 195	CLIH3DeviceNet, 125
GetDataMode	GetDiginValue
CMcsUsbDacqNet, 195	CStg200xBasicNet, 445
GetDDecp	GetDigitalData

CMeaDigitalDataFunctionNet, 285	CStg200xBasicNet, 448
GetDigitalIn	GetEnableAmplifierProtectionSwitch
CPPCFunctionNet, 349	CStg200xBasicNet, 448, 449
CPPS_FunctionNet, 356	GetEnabledChannelsInGroup
GetDigitalOutPortValve	CCMOSMea_FunctionNet, 37
CWarnerValveControllerDeviceNet, 562	CDacqGroupChannelSelectionTemplateNet< Dac
GetDigitalOutPortValveEvent	qGroupChannelEnumTemplateNet, Dac
CWarnerValveControllerDeviceNet, 578	qGroupChannelEnumTemplate, CDevice
GetDigitalPortDirection	GroupChannelInfoTemplateNet >, 54
CWarnerValveControllerDeviceNet, 563	GetEnableHeaterLimit
GetDigitalPortDirectionEvent	CTcxDeviceNet, 502
CWarnerValveControllerDeviceNet, 578	GetEnableThermocouple
GetDigitalSource	CTcxDeviceNet, 502
CMcsUsbDacqNet, 196, 197	GetEntry
GetDigitalStimulatorTrigger	CMcsUsbListEntryNet, 234
CW2100_StimulatorFunctionNet, 534	GetEntryCount
GetDigitalStimulatorTriggerSlope	CMcsUsbListEntryNet, 235
CW2100_StimulatorFunctionNet, 534	GetEnumerationSpeed
GetDigout	CMeaDeviceNet, 280
CFluidControlDeviceNet, 76	GetErrorAirpressure
CRoboDacqNet, 384	CRoboDeviceNet, 396
GetDigoutMode	GetErrorCurrentAirvalve
CStg200xBasicNet, 445	CRoboDeviceNet, 396
GetDigoutValue	GetErrorText
CStg200xBasicNet, 446	CMcsUsbNet, 248
GetDIO	GetErrorVoltage12V
CMcsBus_FYIExtensionNet, 132	CRoboDeviceNet, 396
GetDischargeCapacity	GetErrorVoltage5V
CMultiBatteryChargerDeviceNet, 301	CRoboDeviceNet, 396
GetDischargeCurrent	GetErrorVoltageAirvalve
CMultiBatteryChargerDeviceNet, 301	CRoboDeviceNet, 396
GetDischargeCurrentSetPoint	GetErrorVoltageRs485A
CMultiBatteryChargerDeviceNet, 301	CRoboDeviceNet, 396
GetDisplayMode	GetErrorVoltageRs485B
CWarnerValveControllerDeviceNet, 563	CRoboDeviceNet, 396
GetDisplayModeEvent	GetErrorVoltageValves
CWarnerValveControllerDeviceNet, 578	CRoboDeviceNet, 396
GetDisplayText	GetEventData
CRoboDacqNet, 384	CPositionIIDeviceNet, 335
GetDMax	GetExternalElectrodeEnable
CTcxDeviceNet, 502	CStg200xBasicNet, 449
GetDMin	GetFAAmplification
CTcxDeviceNet, 502	CStg200xBasicNet, 450
GetDownsampleFactor	GetFilter
CRoboDacqNet, 384	CRoboDacqNet, 384
GetDSPHighPassByIndex	GetFilterAttributes
CIntanMea_FunctionNet, 115	CFilterConfigurationNet, 69
GetDuration	GetFilterCoeffs
CMeaCoatDeviceNet, 272	CRoboDacqNet, 384
GetDuty	GetFilterParametersHeadstage
CTcxDeviceNet, 502	CWClassicFunctionNet, 584
GetEEpromPage	GetFilterProperties
CLIH3DeviceNet, 125	CSCUFunctionNet, 425
GetElectrodeDacMux	CW2100_FunctionNet, 527
CStg200xBasicNet, 446	GetFilterProperty
GetElectrodeEnable	· ·
	CMcsUsbDacqNet, 198 CSCUFunctionNet, 425
CStg200xBasicNet, 447	
GetElectrodeMode	CW2100_FunctionNet, 527

O :E: ID: 1	00110011
GetFinalDischargeVoltage	CCMOSMea_FunctionNet, 40
CMultiBatteryChargerDeviceNet, 302	GetGroupID
GetFirmwareVersion	CCMOSMea_FunctionNet, 40
CMcsUsbNet, 248	CDacqGroupChannelSelectionTemplateNet< Dac
GetFirmwareVersionFromFile	qGroupChannelEnumTemplateNet, Dac
CMcsUsbFactoryNet, 225	qGroupChannelEnumTemplate, CDevice
GetFirmwareVersionFromHexFile	GroupChannelInfoTemplateNet >, 54, 55
CMcsUsbFactoryNet, 225	GetGroupNumberOfChannels
GetFPGAFirmwareType	CCMOSMea_FunctionNet, 41
CW2100_FunctionNet, 527	CDacqGroupChannelSelectionTemplateNet< Dac
GetFrameErrorCounter	qGroupChannelEnumTemplateNet, Dac
CTEERFunctionNet, 515	qGroupChannelEnumTemplate, CDevice
GetFrequency	GroupChannelInfoTemplateNet $>$, 55
CRadioControledDevicesNet, 366	GetGroupResolutionPerDigit
GetFrequencyRange	CCMOSMea_FunctionNet, 41
CPgaDeviceNet, 331	GetGroupSampleSize
GetGain	CCMOSMea_FunctionNet, 41
CMeaDeviceNet, 280	CDacqGroupChannelSelectionTemplateNet< Dac
CPgaDeviceNet, 332	qGroupChannelEnumTemplateNet, Dac
GetGate	qGroupChannelEnumTemplate, CDevice
CCMOSMea_FunctionNet, 37	GroupChannelInfoTemplateNet >, 55
GetGilsonDevice	GetGroupType
CRoboocyte2DeviceNet, 410	CCMOSMea_FunctionNet, 41
GetGlobalRepeat	CDacqGroupChannelSelectionTemplateNet< Dac
CDigOutStimulatorFunctionNet, 60	qGroupChannelEnumTemplateNet, Dac
GetGNDI	qGroupChannelEnumTemplate, CDevice
CCMOSMea_FunctionNet, 37	GroupChannelInfoTemplateNet >, 55
GetGroupADCBits	GetGroupUnit
CCMOSMea_FunctionNet, 37	CCMOSMea_FunctionNet, 41, 42
GetGroupChannelBitmaskBySelect	GetGyroRange
CCMOSMea_FunctionNet, 38	CW2100_FunctionNet, 527
GetGroupChannelBitmaskHS1NCBathCurrent	GetHardwareMaxRange
CCMOSMea_FunctionNet, 38	CMcsUsbDacqNet, 200
GetGroupChannelBitmaskHS1NCCol2Current	GetHardwareMinRange
CCMOSMea_FunctionNet, 38	CMcsUsbDacqNet, 200
GetGroupChannelBitmaskHS1NChipTemp	GetHardwareRevision
CCMOSMea_FunctionNet, 38	CMcsUsbNet, 248
GetGroupChannelBitmaskHS1Sidebands	GetHasChecksum
CCMOSMea_FunctionNet, 39	CWClassicFunctionNet, 584
GetGroupChannelBitmaskHS1TriggerStatus	GetHashCode
CCMOSMea_FunctionNet, 39	HeadstageIDTypeObject, 603
GetGroupChannelBitmaskIFDigChannels	GetHasRedLedHeadstage
CCMOSMea_FunctionNet, 39	CWClassicFunctionNet, 584
GetGroupChannelBitmaskInterfaceADC	GetHasThermocouple
CCMOSMea_FunctionNet, 39	CTcxDeviceNet, 502
GetGroupChannelBitmaskPacketFrameContext	GetHeadstage
CCMOSMea_FunctionNet, 40	CStg200xBasicNet, 450
GetGroupChannelBitmaskSTG1DACSignal	GetHeadstageActive
CCMOSMea_FunctionNet, 40	CMcsUsbNet, 249
GetGroupChannelDatal16	GetHeadstageAdcBits
CMcsUsbDacqNet, 198	CSCUFunctionNet, 425
GetGroupChannelDatal32	GetHeadstageAdcRangeInMicroVolt
CMcsUsbDacqNet, 198	CSCUFunctionNet, 426
GetGroupChannelDataUI16	GetHeadstageDacBits
CMcsUsbDacqNet, 199	CSCUFunctionNet, 426
GetGroupChannelDataUl32	GetHeadstageDacCurrentRangeInMicroAmpere
CMcsUsbDacqNet, 199	CSCUFunctionNet, 426
GetGroupDCOffset	GetHeadstageDacCurrentResolutionInNanoAmpere

CSCUFunctionNet, 427	CRoboDacqNet, 384
GetHeadstageDacVoltageRangeInMilliVolt	GetIMax
CSCUFunctionNet, 427	CTcxDeviceNet, 503
GetHeadstageDacVoltageResolutionInMicroVolt	GetIMin
CSCUFunctionNet, 427	CTcxDeviceNet, 503
GetHeadstageGainInPermille	GetImpedanceResult
CSCUFunctionNet, 428	CIntanMea FunctionNet, 115
GetHeadstageID	GetImpedanceTestFrequency
CMcsUsbNet, 249	CMealmpedanceDeviceNet, 291
CSCUFunctionNet, 428	GetImpId
GetHeadstageNumberOfAnalogChannels	CPositionImpDeviceNet, 343
CSCUFunctionNet, 428	GetImplantCurrentSetpoint
GetHeadstageNumberOfStimulationChannels	CPositionIIDeviceNet, 335
CSCUFunctionNet, 429	GetImplantResult
GetHeadstageOnOff	CPositionIIDeviceNet, 336
CW2100_FunctionNet, 528	GetImplantState
CWClassicFunctionNet, 584	CPositionIIDeviceNet, 336
GetHeadstagePresent	GetInMovement
CMcsUsbNet, 249	CRoboDeviceNet, 396
GetHeadstageSamplerate	GetIntanRegister
CSCUFunctionNet, 429	CIntanNea FunctionNet, 115
GetHeadstageSamplingActive	GetIntBuffer
CW2100_FunctionNet, 528	CGenericDevelopDeviceNet, 90
GetHeadstageSerialNumber	GetIO
CSCUFunctionNet, 429	CWarnerValveControllerDeviceTesterFunctionNet,
GetHeaterLimit	581
CTcxDeviceNet, 503	GetlOut 500
GetHeaterTemp	CTcxDeviceNet, 503
CTcxDeviceNet, 503	GetloVoltage
GetHighCurrentRange	CInterfaceboard2FunctionNet, 117
CWarnerUssingFunctionNet, 546	GetLastAnswer
GetHighpassFilterEnable	CGilsonDeviceNet, 101
CFilterConfigurationNet, 69	GetLastUSBError
GetHWConfig	CMcsUsbNet, 250
CRoboDeviceNet::RoboMainLowLevelCommands,	GetLatency
609	CMcsBus_SensorNet, 156
GetHWRevision	GetLatencyCounter
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_SensorNet, 156
609	GetLayoutConfiguration
GetHWRevisionEeprom	CMEA2100x256FunctionNet, 262
CMcsBusNet, 170	GetLEDSwitch
Getl	CMcsBus_ExtensionNet, 131
CTcxDeviceNet, 503	GetLength
GetIC	CRobo_FYIProgram_FunctionNet, 374
CRoboDacqNet, 384	GetLiquidResistance
GetlClamp	CTEERFunctionNet, 515
CRoboDacqNet, 384	CWarnerUssingFunctionNet, 546
GetlCoeff	GetListmodeIndexRange
CRobo_FYITemp_FunctionNet, 376	CStg200xBasicNet, 450
GetICOffset	GetListmodeTriggerSource
CRoboDacqNet, 384	CStg200xBasicNet, 450
GetIDecp	GetLowCurrentRange
CTcxDeviceNet, 503	CWarnerUssingFunctionNet, 547
GetIdent	GetLowerFrequencyByIndex
CMcsUsbNet, 249	CIntanMea_FunctionNet, 115
GetIdleModeOffset	GetMajor
CWarnerUssingFunctionNet, 546	DriverVersionNet, 594
GetlGain	GetMaxChunkSize Byte
GUIGAIII	GOUNGAOHUHKOIZO DYLO

CTEERFunctionNet, 516	GetMCCurrentEeprom
GetMaxCurrent	CMcsBus_MotorControlNet, 138
CMeaCoatDeviceNet, 272	GetMCCurrentMode
GetMaxDurationHighCurrentInMicroSec	CMcsBus_MotorControlNet, 138
CMultiwellOptoStimFunctionNet, 314	GetMCCurrentModeEeprom
GetMaxDutyCycleHighCurrent	CMcsBus MotorControlNet, 138
CMultiwellOptoStimFunctionNet, 314	GetMCCurrentModeShortCommand
GetMaxHeaterPowerMultiwell	CMcsBus_MotorControlNet, 138
CTcxDeviceNet, 504	GetMCCurrentPosition
GetMaxNoPressure	CMcsBus_MotorControlNet, 138
CRoboDeviceNet::RoboMainLowLevelCommands,	GetMCCurrentShortCommand
609	CMcsBus_MotorControlNet, 138
GetMaxNoPressureWaitTime	GetMCCurrentSpeed
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_MotorControlNet, 138
609	GetMCMaxAcceleration
GetMaxNumberOfHeadstages	CMcsBus_MotorControlNet, 139
CSCUFunctionNet, 430	GetMCMaxAccelerationEeprom
GetMaxNumOfColumns	CMcsBus_MotorControlNet, 139
CCMOSMea FunctionNet, 42	GetMCMaxCurrent
GetMaxP	CMcsBus_MotorControlNet, 139
CTcxDeviceNet, 504	GetMCMaxCurrentEeprom
GetMaxpDecp	CMcsBus_MotorControlNet, 139
CTcxDeviceNet, 504	GetMCMaxSpeed
GetMaxpMax	CMcsBus_MotorControlNet, 139
CTcxDeviceNet, 504	GetMCMaxSpeedEeprom
GetMaxpMin	CMcsBus_MotorControlNet, 139
CTcxDeviceNet, 504	GetMCMaxTravel
GetMaxPower	CMcsBus_MotorControlNet, 139
COkuvisionStimulatorDeviceNet, 323	GetMCMaxTravelEeprom
CRobo_FYITemp_FunctionNet, 376	CMcsBus_MotorControlNet, 140
GetMaxPressureWaitTime	GetMCMaxTravelShortCommand
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_MotorControlNet, 140
609	GetMCMovement
GetMaxSamplingFrequency	CMcsBus_MotorControlNet, 140
CMcsUsbDacqNet, 200	GetMCNewPosition
GetMaxStimulusChannelsPerHeadstage	CMcsBus_MotorControlNet, 140
CSCUFunctionNet, 430	GetMCOutputOnOff
GetMaxVoltage	CMcsBus_MotorControlNet, 140
CMeaCleanDeviceNet, 268	GetMCPhase
COkuvisionStimulatorDeviceNet, 323	CMcsBus_MotorControlNet, 140
GetMCAcceleration	GetMCPhaseOffset
CMcsBus MotorControlNet, 136	CMcsBus_MotorControlNet, 140
GetMCAccelerationEeprom	GetMCReference
CMcsBus_MotorControlNet, 136	CMcsBus MotorControlNet, 141
GetMCAccelerationShortCommand	GetMCReferenceCurrent
CMcsBus_MotorControlNet, 137	
	CMcsBus_MotorControlNet, 141
GetMCAxisRevisionEeprom	GetMCReferenceCurrentEeprom
CMcsBus_MotorControlNet, 137	CMcsBus_MotorControlNet, 141
GetMCBreakCurrent CMccBus MeterControlNet 127	GetMCRegulatorGain
CMcSBus_MotorControlNet, 137	CMcsBus_MotorControlNet, 141
GetMCBreakCurrentEeprom	GetMCRegulatorGainEeprom
CMcsBus_MotorControlNet, 137	CMcsBus_MotorControlNet, 141
GetMCConfig CMacRus MaterCentralNet 127	GetMcsBus_Extension
CMcSBus_MotorControlNet, 137	CRoboocyte2DeviceNet, 410
GetMCConfigEeprom	GetMCScalingFactor
CMcsBus_MotorControlNet, 137	CMcsBus_MotorControlNet, 141
GetMCCurrent OMag Para Mater Country IN at 487	GetMCScalingFactorEeprom
CMcsBus_MotorControlNet, 137	CMcsBus_MotorControlNet, 141

GetMCSpeed	GetNeurochipMemoryData
CMcsBus_MotorControlNet, 142	CCMOSMea_FunctionNet, 42
GetMCSpeedEeprom	GetNeurochipMemorySize
CMcsBus_MotorControlNet, 142	CCMOSMea_FunctionNet, 42
GetMCSpeedShortCommand	GetNIC_MS
CMcsBus_MotorControlNet, 142	CRoboDacqNet, 384
GetMCSpeedUnitEeprom	GetNUC_MS
CMcsBus_MotorControlNet, 142	CRoboDacqNet, 384
GetMCStandbyCurrent	GetNumAmplifications
CMcsBus_MotorControlNet, 142	CPgaDeviceNet, 332
GetMCStandbyCurrentEeprom	GetNumber
CMcsBus_MotorControlNet, 142	CMeaSwitchDeviceNet, 293
GetMCStandbyTime	CSw2to64DeviceNet, 497
CMcsBus_MotorControlNet, 142	GetNumberOfAnalogChannels
GetMCStandbyTimeEeprom	CStg200xBasicNet, 450
CMcsBus_MotorControlNet, 143	CStimulusFunctionNet, 489
GetMea21UsbPort	CW2100 StimulatorFunctionNet, 534
CMcsUsbNet, 250	GetNumberOfAudioChannels
GetMeaLayout	CMeaAudioFunctionNet, 265
CMcsUsbDacqNet, 200	GetNumberOfAvailableChambers
GetMemoryUsageDAC	CWarnerUssingFunctionNet, 547
CStg200xDownloadBasicNet, 471	GetNumberOfAvailableSamples
GetMemoryUsageSyncout	CTEERFunctionNet, 516
CStg200xDownloadBasicNet, 471	GetNumberOfChannels
GetMinimalThreshold	CDigOutStimulatorFunctionNet, 61
CMcsBus_SensorNet, 156	GetNumberOfDataBits
GetMinNoPressureWaitTime	CMcsUsbDacqNet, 201
CRoboDeviceNet::RoboMainLowLevelCommands,	GetNumberOfDevices
609	CMcsUsbListNet, 238
GetMinor	GetNumberOfHardwareSlotsForChambers
DriverVersionNet, 594, 595	CWarnerUssingFunctionNet, 547
GetMinPressure	GetNumberOfHWADCChannels
CRoboDeviceNet, 396	CMcsUsbDacqNet::CHWInfo, 113
CRoboDeviceNet::RoboMainLowLevelCommands,	GetNumberOfHWDACPaths
609	CStg200xBasicNet, 450
GetMinPressureWaitTime	GetNumberOfHWDigitalChannels
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsUsbDacqNet::CHWInfo, 113
610	GetNumberOfStimulationElectrodes
GetMinSamplingFrequencyStepsize	CStg200xBasicNet, 450
CMcsUsbDacqNet, 201	GetNumberOfStimulationSourcesPerElectrode
GetMinVoltage	CStg200xBasicNet, 451
CMeaCleanDeviceNet, 268	GetNumberOfSupportedGroups
GetModeSelect	CCMOSMea_FunctionNet, 42
CPulseGeneratorFunctionNet, 363	CDacqGroupChannelSelectionTemplateNet< Dac-
GetModuleCurrent	qGroupChannelEnumTemplateNet, Dac-
CStg200xDownloadNet, 480	qGroupChannelEnumTemplate, CDevice-
GetModuleTemp	GroupChannelInfoTemplateNet >, 55
CStg200xDownloadNet, 480	GetNumberOfSyncoutChannels
GetMovementError	CStg200xBasicNet, 451
CRoboDeviceNet, 396	CW2100_StimulatorFunctionNet, 535
GetMovePump	GetNumberOfTriggerInputs
CMcsBus_SensorNet, 156	CStg200xBasicNet, 451
GetMultiHeadstageMode	CW2100_StimulatorFunctionNet, 535
CW2100_FunctionNet, 528	GetNumConfigurations
GetMultiplexedDataChannelsInBlock	CMcsUsbNet, 250
CStimulusFunctionNet, 488	GetNumControlChannels
GetNanoVoltsPerKelvin	CTcxDeviceNet, 504
CMcsBus_TempSensorNet, 162	GetNumDestinations

CMcsUsbFactoryNet, 225	GetPlateClampStateByHeadstageEvent
GetNumDevices	CMultiwellCallbackFunctionNet, 306
CTcxDeviceNet, 504	GetPlateMux
GetNumEntries	CMultiwellDeviceNet, 308, 309
DriverVersionNet, 595	GetPlateType
GetNumFrequencyRanges	CMultiwellDeviceNet, 309
CPgaDeviceNet, 332	GetPMax
GetNumMeasureChannels	CTcxDeviceNet, 505
CTcxDeviceNet, 504	GetPMin
GetNUV_MS	CTcxDeviceNet, 505
CRoboDacqNet, 385	GetPoti
GetOffsetCurrent	CMcsUsbDacqNet, 201
CMeaCoatDeviceNet, 273	GetPOut
GetOnOff	CTcxDeviceNet, 505
CPositionIIDeviceNet, 336	GetPowerMuxPlate
CTcxDeviceNet, 505	CMultiwellDeviceNet, 309
GetOutputCurrent	GetPowerStrength
CMeaCoatDeviceNet, 273	CPositionIIDeviceNet, 337
GetOutputRate	GetPressure
CStg200xBasicNet, 451	CMcsBus_SensorNet, 156
GetOutputVoltage	GetPressureOffset
CMeaCleanDeviceNet, 268	CMcsBus_SensorNet, 157
GetP	GetPressureRange
CTcxDeviceNet, 505	CPPCFunctionNet, 349
GetParameter	GetPulseform
CRoboDeviceNet::RoboMainLowLevelCommands,	COkuvisionStimulatorDeviceNet, 323
610	GetPulseLength
GetPattern CMacSwitzh PaviacNet 200	CPulseGeneratorFunctionNet, 363
CMeaSwitchDeviceNet, 293	GetPumpCouple
GetPatternBool CMeaSwitchDeviceNet, 294	CPPS_FunctionNet, 356 GetPumpEnableSpeedRatio
GetPauseDuration	CPPS_FunctionNet, 356
CMeaCoatDeviceNet, 273	GetPumpFastOnOff
GetPCoeff	CPPS_FunctionNet, 356
CRobo_FYITemp_FunctionNet, 376	GetPumpFastSpeed
GetPDecp	CPPS FunctionNet, 356
CTcxDeviceNet, 505	GetPumpFunctionSpeeds
GetPeriod	CPPS FunctionNet, 357
CPulseGeneratorFunctionNet, 363	GetPumpManualOnOff
GetPeriod us	CPPS FunctionNet, 357
CTEERFunctionNet, 516	GetPumpMaxSpeed
GetPermanentCurrentInMicroAmp	CPPS_FunctionNet, 357
CMultiwellOptoStimFunctionNet, 315	GetPumpModeType
GetPGain	CPPCFunctionNet, 349
CRoboDacqNet, 385	CPPS_FunctionNet, 357
GetPhases	GetPumpSpeed
CRoboDeviceNet::RoboMainLowLevelCommands,	CRoboFluidDeviceNet, 407
610	GetPumpSpeedRatio
GetPicFirmwareType	CPPS FunctionNet, 357
CW2100_FunctionNet, 528	GetPumpSpeedUnit
GetPiezoState	CPPCFunctionNet, 350
CMcsBus_SensorNet, 156	CPPS_FunctionNet, 357
GetPlateClampLockState	GetPWM
CMultiwellDeviceNet, 308	CFluidControlDeviceNet, 76
GetPlateClampState	GetPwrOut
CMultiwellDeviceNet, 308	CTcxDeviceNet, 505
GetPlateClampStateByHeadstage	GetPwrSet
CMultiwellCallbackFunctionNet, 305	CTcxDeviceNet, 506

GetRatedCapacity	GetRotatePump
CMultiBatteryChargerDeviceNet, 302	CMcsBus_SensorNet, 157
GetReady	GetROut
CMealmpedanceDeviceNet, 291	CTcxDeviceNet, 506
GetRecordingNumber	GetRTC
CRoboDacqNet, 385	COkuvisionStimulatorDeviceNet, 324
GetReferenceElectrodeMode	CPositionIIDeviceNet, 337
CSCUFunctionNet, 430	GetSampleBufferChunk
GetReferenceElectrodeSwitchState	CTEERFunctionNet, 516
CSCUFunctionNet, 430	GetSampleInterval
GetReferenceTemperature	CLIH3DeviceNet, 125
CFluidControlDeviceNet, 76	GetSamplePeriode
GetRegulationTimeouts	CMcsBus_SensorNet, 158
CMcsBus_SensorNet, 157	GetSampleRate
GetRegulatorFactor	CTEERFunctionNet, 517
CMcsBus_SensorNet, 157	GetSamplerate
	•
GetRegulatorOnOff	CMcsUsbDacqNet, 201
CMcsBus_SensorNet, 157	GetSampleVoltageBuffer_uV
CRobo_FYITemp_FunctionNet, 376	CTEERFunctionNet, 517
GetRegulatorStatus	GetScaleFactorU1
CMcsBus_SensorNet, 157	CTEERFunctionNet, 517
GetRepeats	GetScaleFactorU2
CProgramPressureCurveNet, 361	CTEERFunctionNet, 517
GetRes1	GetScanHeadstagesResult
CTcxDeviceNet, 506	CWClassicFunctionNet, 585
GetRes2	GetScreen
CTcxDeviceNet, 506	CRoboDacqNet, 385
GetResetFilter	GetSearchReferenceFastAccel
CWClassicFunctionNet, 585	CRoboDeviceNet::RoboMainLowLevelCommands,
GetResistanceC	610
CRoboDacqNet, 385	GetSearchReferenceFastSpeed
GetResistanceV	CRoboDeviceNet::RoboMainLowLevelCommands,
CRoboDacqNet, 385	610
GetResolutionPerDigit	GetSearchReferenceFineAccel
CMcsUsbDacqNet, 201	CRoboDeviceNet::RoboMainLowLevelCommands,
GetResS	610
CTcxDeviceNet, 506	GetSearchReferenceFineSpeed
GetResult	CRoboDeviceNet::RoboMainLowLevelCommands,
CMealmpedanceDeviceNet, 291	610
GetResX	GetSearchReferenceMethod
CTcxDeviceNet, 506	CRoboDeviceNet::RoboMainLowLevelCommands,
GetRFConnectionStatus	611
CWClassicFunctionNet, 585	GetSearchReferenceMoveOut
GetRFFrequency	CRoboDeviceNet::RoboMainLowLevelCommands,
CPositionImpDeviceNet, 343	611
GetRFFrequencyHeadstage	GetSearchReferenceOffsetPos
CWClassicFunctionNet, 585	CRoboDeviceNet::RoboMainLowLevelCommands.
GetRFFrequencyReceiver	611
CWClassicFunctionNet, 585	GetSelectedChannels
GetRFPower	CW2100_FunctionNet, 528
CWClassicFunctionNet, 585	GetSelectedHeadstage
GetRoboDacq	CWClassicFunctionNet, 585
CRoboocyte2DeviceNet, 410	GetSelectedHeadstageState
GetRoboFluidDevice	CW2100_FunctionNet, 528
CEncapsulatorDeviceNet, 64	GetSensorType
CRoboocyte2DeviceNet, 410	CTcxDeviceNet, 506
GetRotaryPositionCode	GetSerialNumber
CTEERFunctionNet, 516	CMcsUsbNet, 250
C. ELIT GIOGOTITO, CTO	C.11000001101, 200

GetSerialNumberHeadstage	GetStimulusParametersCache
CWClassicFunctionNet, 585	CW2100_FunctionNet, 528
GetSetpoint	GetStimulusParametersFromSelectedHS
CTcxDeviceNet, 506	CW2100_FunctionNet, 528
GetSetpointDecp	GetStimulusSites
CTcxDeviceNet, 507	CCMOSMea_FunctionNet, 43
GetSetpointMax	GetStiumlusParameters
CTcxDeviceNet, 507	CW2100_FunctionNet, 529
GetSetpointMin	GetStopTriggerSlope
CTcxDeviceNet, 507	CDigOutStimulatorFunctionNet, 61
GetShortBuffer	GetSubChannel
CGenericDevelopDeviceNet, 91	CMcsBus_MotorControlNet, 143
GetSimulation	GetSupplyVoltage
CRoboDacqNet, 385	CPPCFunctionNet, 350
GetSingleHeater	CPPS_FunctionNet, 357
CMcsBus_FYIExtensionNet, 132	GetSweepCount
GetSingleValve	CStg200xDownloadBasicNet, 472
CFluidControlDeviceNet, 76	GetSync
CRoboFluidDeviceNet, 407	CWarnerValveControllerDeviceTesterFunctionNet
GetSlope	581
CMeaCleanDeviceNet, 269	GetSyncoutMap
CMeaCoatDeviceNet, 273	CStg200xBasicNet, 452
GetSoftwareKey	GetSyncState
CMcsUsbNet, 250	CMcsBus_SensorNet, 158
GetSoftwareKeyString	GetTableName
CMcsUsbNet, 250	CWarnerValveControllerDeviceNet, 563
GetSollPressure	GetTableNamebyIndex
CMcsBus_SensorNet, 158	CWarnerValveControllerDeviceNet, 563
GetSollTemp	GetTableNamebyIndexEvent
CRobo_FYITemp_FunctionNet, 376	CWarnerValveControllerDeviceNet, 578
GetSourceBulk	GetTablepointer
CCMOSMea_FunctionNet, 42	CRetinaLedDeviceNet, 368
GetSourceDrain	GetTemperatur
CCMOSMea FunctionNet, 42	CMcsBus_TempSensorNet, 162
GetSourceGate	GetTestMode
CCMOSMea FunctionNet, 42	CRFFunctionNet, 372
GetStartTriggerSlope	GetThermocoupleCalibration
CDigOutStimulatorFunctionNet, 61	CFluidControlDeviceNet, 77
GetState	CTcxDeviceNet, 507
CRFFunctionNet, 372	GetThermocoupleNanovoltPerKelvin
CRobo FYIProgram FunctionNet, 375	CFluidControlDeviceNet, 77
GetStateDebugData	CTcxDeviceNet, 507
CPositionIIDeviceNet, 337	GetThermocoupleReferenceTemp
GetStateEventData	CTcxDeviceNet, 507
CPositionIIDeviceNet, 338	GetThermocoupleTemp
GetStatus	CTcxDeviceNet, 507
CMcsUsbNet, 250	GetThermocoupleTempAbs
DriverVersionNet, 595	CTcxDeviceNet, 508
GetStatusOfLastCommand	GetThermocoupleTemperature
CMcsUsbNet, 251	CFluidControlDeviceNet, 77
GetStgProgramInfo	GetThermoOffset
CStg200xBasicNet, 451, 452	CMcsBus_TempSensorNet, 162
GetStgVersionInfo	GetThermoTemp
CStg200xBasicNet, 452	CMcsBus_TempSensorNet, 162
GetStimulationPatternMemory	GetThermoVoltage
CW2100_StimulatorFunctionNet, 535	CMcsBus_TempSensorNet, 162
GetStimulatorStatus	GetTimeInPause
COkuvisionStimulatorDeviceNet. 324	CMeaCoatDeviceNet. 273
JONAVISIONISINNALUNDEVICENTI. JAH	OIVICACOALDEVICEINEL. 4/0

GetTimeInPlateau	CIntanMea FunctionNet, 116
CMeaCoatDeviceNet, 274	GetUptimeSeconds
GetTimeResolutionInNanoSeconds	CTEERFunctionNet, 518
CW2100_StimulatorFunctionNet, 535	CWarnerUssingFunctionNet, 550
GetTimeSlot	GetUSBDeviceIDFromFX3Image
CW2100 StimulatorFunctionNet, 535	CMcsUsbFactoryNet, 226
GetTotalMemory	GetUsbListEntries
•	
CStg200xBasicNet, 453 CStimulusFunctionNet, 489	CMcsUsbListNet, 238
•	GetUsbListEntry
GetTotalNumberOfDigitalPorts	CMcsUsbListNet, 238
CWarnerValveControllerDeviceNet, 564	CMcsUsbNet, 251
GetTotalNumberOfTables	GetUseBubble
CWarnerValveControllerDeviceNet, 564	CPPS_FunctionNet, 357
GetTotalNumberOfValves	GetUsercodeFromBitFile
CWarnerValveControllerDeviceNet, 564	CMcsUsbFactoryNet, 226
GetTotalTableSize	GetUsercodeFromFlash
CWarnerValveControllerDeviceNet, 564	CMcsUsbFactoryNet, 226
GetTrigger	GetUserDefinedName
CStg200xDownloadBasicNet, 472	CW2100_FunctionNet, 529
GetTriggerSource	GetUserDefinedNameCache
CStg200xBasicNet, 453	CW2100_FunctionNet, 529
GetU1Offset	GetUserDefinedNameFromSelectedHS
CWarnerUssingFunctionNet, 547	CW2100_FunctionNet, 529
GetU1Reference	GetUserParameter
CWarnerUssingFunctionNet, 548	CRoboDeviceNet::RoboMainLowLevelCommands,
GetU2Offset	611
CWarnerUssingFunctionNet, 548	GetUShortBuffer
GetU2Reference	CGenericDevelopDeviceNet, 94
CWarnerUssingFunctionNet, 548	GetUV
GetUByteBuffer	CRoboDacqNet, 386
CGenericDevelopDeviceNet, 91	GetUVOffset
GetUC	CRoboDacqNet, 386
CRoboDacqNet, 385	GetValue
GetUClamp	CTcxDeviceNet, 508
CRoboDacqNet, 385	GetValueHires
GetUCOffset	CTcxDeviceNet, 508
CRoboDacqNet, 385	GetValve
GetUintA	CFluidControlDeviceNet, 78
CFilterCoefficientsNet, 67	CRoboFluidDeviceNet, 407
GetUintB	GetValve1
CFilterCoefficientsNet, 67	CRobo FYIProgram FunctionNet, 375
GetUIntBuffer	GetValve2
CGenericDevelopDeviceNet, 93	
GetUnit	CRobo_FYIProgram_FunctionNet, 375
	GetValveActive
CTcxDeviceNet, 508	CPPCFunctionNet, 350
GetUnitDescription	CWarnerValveControllerDeviceNet, 564
CWarnerUssingFunctionNet, 549	GetValveActiveEvent
GetUnitExponent	CWarnerValveControllerDeviceNet, 578
CWarnerUssingFunctionNet, 549	GetValveBoardRevision
GetUnitName	CWarnerValveControllerDeviceNet, 565
CWarnerUssingFunctionNet, 549	GetValveBoardRevisionEvent
GetUnitsPerDigit	CWarnerValveControllerDeviceNet, 578
CWarnerUssingFunctionNet, 550	GetValveBoardRevisionString
GetUOut	CWarnerValveControllerDeviceNet, 565
CTcxDeviceNet, 508	GetValveDigitalInPort
GetUpdateDisplay	CWarnerValveControllerDeviceNet, 565
CRoboDacqNet, 385	GetValveDigitaIInPortEvent
GetUpperFrequencyByIndex	CWarnerValveControllerDeviceNet, 578

GetValveLedOn	GetVoltage
CWarnerValveControllerDeviceNet, 565	COkuvisionStimulatorDeviceNet, 324
GetValveLedOnEvent	GetVoltage12V
CWarnerValveControllerDeviceNet, 579	CRoboDeviceNet, 397
GetValveManualGroup	GetVoltage12VLimit
CWarnerValveControllerDeviceNet, 566	CRoboDeviceNet, 397
GetValveManualGroupEvent	GetVoltage5V
CWarnerValveControllerDeviceNet, 579	CRoboDeviceNet, 397
GetValveManualState	GetVoltage5VLimit
CWarnerValveControllerDeviceNet, 566	CRoboDeviceNet, 397
GetValveManualStateEvent	GetVoltageAirvalve
CWarnerValveControllerDeviceNet, 579	CRoboDeviceNet, 397
GetValveMode	GetVoltageAirvalveLimit
CWarnerValveControllerDeviceNet, 566	CRoboDeviceNet, 397
GetValveModeEvent	GetVoltageClampControllerParam_D
CWarnerValveControllerDeviceNet, 579	CWarnerUssingFunctionNet, 550
GetValves	GetVoltageClampControllerParam_I
CMcsBus_FYIExtensionNet, 132	CWarnerUssingFunctionNet, 551
GetValvesActiveMap	GetVoltageClampControllerParam_P
CWarnerValveControllerDeviceNet, 566	CWarnerUssingFunctionNet, 551
GetValvesManualStateMap	GetVoltageRange
CWarnerValveControllerDeviceNet, 567	CGrapheneFunctionNet, 105
GetValveTableEntry	GetVoltageRangeIndex
CWarnerValveControllerDeviceNet, 567	CMcsUsbDacqNet, 201
GetVDD3I	GetVoltageRangeInMicroVolt
CCMOSMea_FunctionNet, 43	CMcsUsbDacqNet, 202
GetVDDI	CStg200xBasicNet, 453
CCMOSMea_FunctionNet, 43	CStimulusFunctionNet, 489
GetVdVs	CW2100_StimulatorFunctionNet, 535
CGrapheneFunctionNet, 104	GetVoltageRangeInMilliVolt
GetVdVsDAC	CMcsUsbDacqNet, 202
CGrapheneFunctionNet, 104, 105	GetVoltageReached
GetVersion	CGrapheneFunctionNet, 105
CMcsUsbNet, 251	GetVoltageResolution
GetVersionInt	CGrapheneFunctionNet, 106
DriverVersionNet, 595, 596	GetVoltageResolutionInMicroVolt
GetVersionString	CStg200xBasicNet, 453
DriverVersionNet, 596	CStimulusFunctionNet, 489
GetVMMaxNegativeCurrent	CW2100_StimulatorFunctionNet, 535
CMcsBus_VoltageModeNet, 164	GetVoltageRs485A
GetVMMaxNegativeCurrentEeprom	CRoboDeviceNet, 397
CMcsBus_VoltageModeNet, 164	GetVoltageRs485ALimit
GetVMMaxNegativeVoltage	CRoboDeviceNet, 397
CMcsBus_VoltageModeNet, 164	GetVoltageRs485B
GetVMMaxNegativeVoltageEeprom	CRoboDeviceNet, 397
CMcsBus_VoltageModeNet, 164	GetVoltageRs485BLimit
GetVMMaxPositiveCurrent	CRoboDeviceNet, 398
CMcsBus_VoltageModeNet, 165	GetVoltageValves
GetVMMaxPositiveCurrentEeprom	CRoboDeviceNet, 398
CMcsBus_VoltageModeNet, 165	GetVoltageValvesLimit
GetVMMaxPositiveVoltage	CRoboDeviceNet, 398
CMcsBus_VoltageModeNet, 165	GetVolti
GetVMMaxPositiveVoltageEeprom	CTcxDeviceNet, 508
CMcsBus_VoltageModeNet, 165	GetWaveform
GetVMOutputOnOff	CTEERFunctionNet, 518
CMcsBus_VoltageModeNet, 165	GetWaveLengthInNanometer
GetVMVoltage	CMultiwellOptoStimFunctionNet, 315
CMcsBus_VoltageModeNet, 165	GetWorkingFrequency

CRFFunctionNet, 373	StimulusParameters, 602
GetWPADebugMode	ToString, 601
CWClassicFunctionNet, 585	Type, 602
GetWPAType	TypeValue, 602
CWClassicFunctionNet, 585	Unknown, 601
GetXGain	UserDefinedName, 602
CRoboDacqNet, 386	Valid, 602
GetXilinxFlashOffset	W16IsW14, 602
CMcsUsbFactoryNet, 226	HeadstageIDTypeObject, 603
GetXilinxFlashReadCommand	AdditionalText, 603
CMcsUsbFactoryNet, 226	_Additional Text, 000 _IdType, 604
GND_SWITCH_BIT	AdditionalText, 604
CW2100_StimulatorFunctionNet, 539	Equals, 603
GroupID	GetHashCode, 603
CDeviceGroupChannelInfoTemplateNet< Dacq-	HeadstageIDTypeObject, 603
GroupChannelEnumTemplateNet >, 58	IdType, 604
GroupType	ToString, 603
CDeviceGroupChannelInfoTemplateNet< Dacq-	HeadStageIDTypeState, 604
GroupChannelEnumTemplateNet $>$, 58	ControlState, 604
	DataState, 604
HasAnalogOut	IdType, 604
CSCUFunctionNet, 431	State, 604
HasGalvanicIsolation	HeadstageType
CSCUFunctionNet, 431	HeadStageIDType, 602
HasHSPowerSwitch	HeadstageTypeEnum
CSCUFunctionNet, 431	HeadStageIDType, 600
HasIMU	HLADacq
HeadStageIDType, 601	CHLADeviceNet, 112
HasOptoCurrentMessurement	HWInfo
HeadStageIDType, 601	
HasRadioControl	CMcsUsbDacqNet, 202
CRadioControledDevicesNet, 366	HwVersion
HasRef	CMcsUsbListEntryNet, 236
CRoboDeviceNet::RoboMainLowLevelCommands,	ID
612	HeadStageIDType, 602
HasRefl	IdProduct
CRoboStatorDeviceNet, 413	DeviceIdNet, 590
HasRefXY	IdType
CRoboStatorDeviceNet, 413	HeadstageIDTypeObject, 604
HasRefZ	HeadStageIDTypeState, 604
CRoboStatorDeviceNet, 413	IdVendor
HasSoftwareKey	DeviceIdNet, 590
CMcsUsbNet, 251	IsAnalogOutEnabled
HeadStageIDType, 600	CSCUFunctionNet, 431
CompareTo, 601	IsAutomaticAnalogOut
ElectricalStimulation, 601	CSCUFunctionNet, 432
Entry, 601	IsBusy
Equals, 601	CPPCFunctionNet, 350
HasIMU, 601	IsChamberAvailable
HasOptoCurrentMessurement, 601	CWarnerUssingFunctionNet, 551
HeadStageIDType, 601	IsChipPowered
HeadstageType, 602	•
- ··	CCMOSMea_FunctionNet, 43
HeadstageTypeEnum, 600	IsConnected
ID, 602	CMcsUsbNet, 251
MeasuringOnly, 601	IsDeviceHighSpeed
NumberOfAnalogChannels, 602	CMcsUsbNet, 251
NumberOfStimulationChannels, 602	IsDeviceHighSpeedCapable
OpticalStimulation, 601	CMcsUsbNet, 252
SN, 602	IsDeviceTypeOf

CMcsUsbListNet, 239	CStg200xBasicNet, 454
IsDigitalChannelDedicated	ListOfChangedTriggers
CMcsUsbDacqNet::CHWInfo, 114	StgStatusNet, 618
IsDigitalOutPortInverted	LoadPressure
CWarnerValveControllerDeviceNet, 567	CPPCFunctionNet, 352
IsDigitalOutPortInvertedEvent	LoadUserFirmware
CWarnerValveControllerDeviceNet, 579	CMcsUsbFactoryNet, 226, 227
IsEqual	LoadValveTable
CFilterCoefficientsNet, 67	CWarnerValveControllerDeviceNet, 569
IsExceptionsEnabled	LockPlateClamp
CMcsUsbNet, 252	CMultiwellDeviceNet, 310
IsGateFloating	,
CCMOSMea_FunctionNet, 43	m_Bottom
IsHeadstageAvailable	CCMOSMeaDeviceNet:: CRegionOfInterestRect,
CSCUFunctionNet, 432	367
IsHeadstageAvailableEvent	m_Left
CSCUFunctionNet, 436	CCMOSMeaDeviceNet::CRegionOfInterestRect,
IsHighCurrentMode	367
CWarnerUssingFunctionNet, 552	m_pGilsonDevice
IsHSPowered	CGilsonDeviceNet, 102
CSCUFunctionNet, 432	m_pMcsBus_MotorControlNet
IsInDacqLegacyMode	CRoboFluidDeviceNet, 408
CSCUFunctionNet, 432	m_pMcsUsb
IsInternalCalibrationFinished	CMcsUsbFunctionNet, 231
CTEERFunctionNet, 518	m_pMcsUsbFunction
CWarnerUssingFunctionNet, 552	CMcsUsbFunctionNet, 231
IsPlateTypeValid	m_pRoboFluidDevice
CMultiwellDeviceNet, 310	CRoboFluidDeviceNet, 408
IsPulseEnabled	m_Right
CWarnerUssingFunctionNet, 552	CCMOSMeaDeviceNet::CRegionOfInterestRect,
IsPumpMotorOn	367
CRoboFluidDeviceNet, 407	m_Top
IsRunning	CCMOSMeaDeviceNet::CRegionOfInterestRect,
CMeaCleanDeviceNet, 269	367
CMeaCoatDeviceNet, 274	Manufacturer
IsSamplingFinished	CMcsUsbListEntryNet, 236
CTEERFunctionNet, 518	MaxBitNumber 500
IsUserTriggerEnabled	DigitalSource< digitalsourceenum >, 590
CLIH3DeviceNet, 125	Mcs, 22
IsValveDigitalInInverted	Mcs::Usb, 22
CWarnerValveControllerDeviceNet, 567	DEVICE_NOT_FOUND, 27
IsValveDigitalInInvertedEvent	enCMosMeaChipType, 27
CWarnerValveControllerDeviceNet, 579	EnSTG200x_STATUS, 27
IsValveOpen	nMos16LV, 27
CWarnerValveControllerDeviceNet, 568	nMos32LV, 27 nMos36LN, 27
IsValveOpenEvent	nMos64LN, 27
CWarnerValveControllerDeviceNet, 579	NOT_CONNECTED, 27
IsValveOpenInAnalogMode	
CWarnerValveControllerDeviceNet, 568	OK, 27 OnChannelData, 27
IsValveOpenInAnalogModeEvent	OnDeviceArrivalRemoval, 27
CWarnerValveControllerDeviceNet, 579	OnError, 27
IsValveOpenInDigitalMode	OnMcsUsbDeviceState, 28
CWarnerValveControllerDeviceNet, 568	OnMcsUsbDeviceStateCallback, 28
IsValveOpenInDigitalModeEvent	OnMwPollStatus, 28
CWarnerValveControllerDeviceNet, 580	OnStg200xDataHandler, 28
ListModeSendStart	OnStg200xErrorHandler, 28
CStg200xBasicNet, 454	OnStgPollStatus, 28
ListModeSendStop	OnUpdateFirmwareProgress, 28
	5 opadio: :::::::::::::::::::::::::::::::::::

OnUpdateFirmwareStatusChange, 28	CMeaDeviceNet, 283
RoboStatusEventDelegate, 29	Measure
unknown, 27	CPathIdentDeviceNet, 327
McsBus	MeasureReservoir
CPPCDeviceNet, 345	CPPCFunctionNet, 352
CPPS_DeviceNet, 354	MeasuringOnly
CRoboDeviceNet, 405	HeadStageIDType, 601
MCSBUS1	MeFunctionNet
FirmwareDestinationNames, 598	CMeaDeviceNet, 284
MCSBUS10	mkfilter
FirmwareDestinationNames, 598	mkfilterNet, 605
MCSBUS11	mkfilter_coef_in_one_set
FirmwareDestinationNames, 598	mkfilterNet, 605
MCSBUS12	mkfilter_highpass_coeff
FirmwareDestinationNames, 598	mkfilterNet, 606
MCSBUS13	mkfilter_highpass_frequency_from_coeff
FirmwareDestinationNames, 598	mkfilterNet, 606
MCSBUS2	mkfilter highpass frequency from k
FirmwareDestinationNames, 598	mkfilterNet, 606
MCSBUS3	mkfilter_highpass_k
FirmwareDestinationNames, 598	mkfilterNet, 606
MCSBUS4	mkfilter_MCS
FirmwareDestinationNames, 599	mkfilterNet, 606
MCSBUS5	mkfilter_MCS_k
FirmwareDestinationNames, 599	mkfilterNet, 606, 607
MCSBUS6	mkfilter_normalize_coeffs_int
FirmwareDestinationNames, 599	mkfilterNet, 607
MCSBUS7	mkfilter_normalize_coeffs_short
FirmwareDestinationNames, 599	mkfilterNet, 607
MCSBUS8	mkfilter_normalize_scale_coeffs_int
FirmwareDestinationNames, 599	mkfilterNet, 607
MCSBUS9	mkfilter_scale_coef_in_one_set
FirmwareDestinationNames, 599	mkfilterNet, 607
McsBus_MotorControl	mkfilterNet, 605
CPeristalticPumpDeviceNet, 330	mkfilter, 605
CPPCDeviceNet, 345	mkfilter_coef_in_one_set, 605
CPPS_DeviceNet, 354	mkfilter_highpass_coeff, 606
CRoboDeviceNet, 405	mkfilter_highpass_frequency_from_coeff, 606
CRoboFluidDeviceNet, 408	mkfilter_highpass_frequency_from_k, 606
McsBus Sensor	mkfilter_highpass_k, 606
CPPCDeviceNet, 345	mkfilter_MCS, 606
CPPS DeviceNet, 354	mkfilter_MCS_k, 606, 607
McsBus_VoltageMode	mkfilter normalize coeffs int, 607
CFluidControlDeviceNet, 80	mkfilter_normalize_coeffs_short, 607
	mkfilter normalize scale coeffs int, 607
McsBus_XY	mkfilter_scale_coef_in_one_set, 607
CRoboDeviceNet, 402	MoveAbs
McsBus_ZI	
CRoboDeviceNet, 402	CRoboDeviceNet, 398
McsUsbDeviceStateEvent	MoveAbsl
CMcsUsbDeviceStatePushFunctionNet, 220	CRoboStatorDeviceNet, 413, 414
CMcsUsbDeviceStatePushNet, 221	MoveAbsXY
MCU1	CRoboStatorDeviceNet, 414
FirmwareDestinationNames, 599	MoveAbsZ
MeaAudioFunctionNet	CRoboStatorDeviceNet, 414
CMeaDeviceNet, 283	MultibootGetCypressImageId
MeaDigitalDataFunctionNet	CMcsUsbNet, 252
CMeaDeviceNet, 283	MultibootGetImageId
MeaFeedbackFunctionNet	CMcsUsbNet, 252

MultibootGetSelectedImage	OnGetValveBoardRevision
CMcsUsbNet, 252	CWarnerValveControllerDeviceNet, 570
MultibootSelectImage	OnGetValveDigitalInPort
CMcsUsbNet, 252	CWarnerValveControllerDeviceNet, 570
MwPollStatusEvent	OnGetValveLedOn
CStg200xDownloadNet, 483	CWarnerValveControllerDeviceNet, 570
Odig=00xDownloadivot, 100	OnGetValveManualGroup
nMos16LV	•
Mcs::Usb, 27	CWarnerValveControllerDeviceNet, 570
nMos32LV	OnGetValveManualState
	CWarnerValveControllerDeviceNet, 570
Mcs::Usb, 27	OnGetValveMode
nMos36LN	CWarnerValveControllerDeviceNet, 570
Mcs::Usb, 27	OnIsDigitalOutPortInverted
nMos64LN	CWarnerValveControllerDeviceNet, 571
Mcs::Usb, 27	OnIsHeadstageAvailable
NOT_CONNECTED	CSCUFunctionNet, 433
Mcs::Usb, 27	OnlsValveDigitalInInverted
NumberOfAnalogChannels	CWarnerValveControllerDeviceNet, 571
HeadStageIDType, 602	Onls Valve Open
NumberOfChannels	
	CWarnerValveControllerDeviceNet, 571
·	OnlsValveOpenInAnalogMode
GroupChannelEnumTemplateNet >, 58	CWarnerValveControllerDeviceNet, 571
NumberOfStimulationChannels	OnlsValveOpenInDigitalMode
HeadStageIDType, 602	CWarnerValveControllerDeviceNet, 571
NumCoefSets	OnMcsUsbDeviceState
CCreateFilterNet, 49	Mcs::Usb, 28
	OnMcsUsbDeviceStateCallback
OK	Mcs::Usb, 28
Mcs::Usb, 27	OnMwPollStatus
OnChannelData	Mcs::Usb, 28
Mcs::Usb, 27	OnStg200xDataHandler
OnDeviceArrivalRemoval	Mcs::Usb, 28
Mcs::Usb, 27	
OnError	OnStg200xErrorHandler
Mcs::Usb, 27	Mcs::Usb, 28
OnGetActiveRunningTableNumber	OnStgPollStatus
CWarnerValveControllerDeviceNet, 569	Mcs::Usb, 28
	OnTableEntryChanged
OnGetAnalogThresholdHigh	CWarnerValveControllerDeviceNet, 571
CWarnerValveControllerDeviceNet, 569	OnUpdateFirmwareProgress
OnGetAnalogThresholdLow	Mcs::Usb, 28
CWarnerValveControllerDeviceNet, 569	OnUpdateFirmwareStatusChange
OnGetAnalogVoltage	Mcs::Usb, 28
CWarnerValveControllerDeviceNet, 569	OpenPipe
OnGetAvailableHeadstages	CGenericDevelopDeviceNet, 94
CSCUFunctionNet, 433	OpenPlateClamp
OnGetCurrentNumberOfValves	CMultiwellDeviceNet, 310
CWarnerValveControllerDeviceNet, 569	
OnGetDigitalOutPortValve	operator=
CWarnerValveControllerDeviceNet, 569	DeviceIdNet, 589
OnGetDigitalPortDirection	OpticalStimulation
_	HeadStageIDType, 601
CWarnerValveControllerDeviceNet, 569	Order
OnGetDisplayMode	CCreateFilterNet, 49
CWarnerValveControllerDeviceNet, 570	CFilterPropertyNet, 73
OnGetPlateClampStateByHeadstage	
CMultiwellCallbackFunctionNet, 305	PatternListStart
OnGetTableNamebyIndex	COctoPotDeviceNet, 319
CWarnerValveControllerDeviceNet, 570	PIC
OnGetValveActive	FirmwareDestinationNames, 599
CWarnerValveControllerDeviceNet 570	PIC2

FirmwareDestinationNames, 599	ReadBlockFromIFBGlobalEEprom
PIC3	CMcsUsbFactoryNet, 227
FirmwareDestinationNames, 599	ReadBlockFromNVMEM
PIC4	CMcsUsbFactoryNet, 227
FirmwareDestinationNames, 599	ReadClipping
PollStatusEvent	CLIH3DeviceNet, 126
CStimulusFunctionNet, 495	ReadEepromRegisterPreconfig
CW2100_StimulatorFunctionNet, 539	CMcsUsbNet, 253
PostCommaA	ReadPipe
CFilterCoefficientsNet::s_FilterAttributesNet, 616	CGenericDevelopDeviceNet, 95
PostCommaB	ReadRegister
CFilterCoefficientsNet::s FilterAttributesNet, 616	CMcsUsbNet, 253
PowerChip	ReadRegister32
CCMOSMea FunctionNet, 43	CMcsUsbNet, 253
PowerHS	
	ReadRegisterTimeSlot
CSCUFunctionNet, 433	CMcsUsbNet, 253
PPCFunction CDPCP in No. 15	Receive
CPPCDeviceNet, 345	CSerialPortNet, 437
PPS_Function	ReceiveString
CPPS_DeviceNet, 354	CSerialPortNet, 437
PreCommaA	RemoveSoftwareKey
CFilterCoefficientsNet::s_FilterAttributesNet, 616	CMcsUsbNet, 254
PreCommaB	RescanHeadstage
CFilterCoefficientsNet::s_FilterAttributesNet, 616	CMcsUsbNet, 254
PrepareAndAppendData	ResetAdcOffset
CStg200xDownloadNet, 480	COctoPotDeviceNet, 320
CStimulusFunctionNet, 490	ResetChannelmap
PrepareAndSendData	CWClassicFunctionNet, 586
CStg200xDownloadNet, 481	ResetDacOffset
CStimulusFunctionNet, 491	COctoPotDeviceNet, 320
PrepareChannelData	ResetHighpassFilter
CDigOutStimulatorFunctionNet, 62	CFilterConfigurationNet, 69
PrepareData	ResetPipe
CStimulusFunctionNet, 491	CGenericDevelopDeviceNet, 95
CW2100 StimulatorFunctionNet, 537	ResetStatus
-	
PrepareDataSync	CStg200xDownloadBasicNet, 472
CW2100_StimulatorFunctionNet, 537	RFFunction
Product	CPositionIIDeviceNet, 341
CMcsUsbListEntryNet, 236	RoboDacq
Program	CHiClampDeviceNet, 110
CProgramPressureCurveNet, 361	RoboDevice
PulseGenerator	CSafeISDeviceNet, 420
CW2100_FunctionNet, 531	RoboError_AnotherMaster
PumpOff	CRoboDeviceNet, 402
CRoboFluidDeviceNet, 407	RoboError_Base
PumpOn	CRoboDeviceNet, 402
CRoboFluidDeviceNet, 407	RoboError_CannotEscapeEndSwitch
	CRoboDeviceNet, 402
QueryTriggerstatus	RoboError_CommandAlreadyInProgress
CStg200xDownloadNet, 482	CRoboDeviceNet, 402
,	RoboError_CommandNotPossible
RampStart	CRoboDeviceNet, 402
COctoPotDeviceNet, 320	RoboError_CommunicationTimeout
Read	
CExternDTesterDeviceNet, 65	CRoboDeviceNet, 402
Read2	RoboError_DacqNotReady
CExternDTesterDeviceNet, 65	CRoboDeviceNet, 402
ReadBlockFromFlash	RoboError_DLLMovementTimeout
CMcsUsbFactoryNet, 227	CRoboDeviceNet, 403
OIVIGOODI ACIDI YINGI, 221	

RoboError_FindReferenceMethod	ScanForHeadstages
CRoboDeviceNet, 403	CWClassicFunctionNet, 586
RoboError_GilsonCommandPending	SelectHeadstage
CRoboDeviceNet, 403	CW2100_FunctionNet, 529
RoboError_GilsonTimeout	SelectTimeSlot
CRoboDeviceNet, 403	CW2100_StimulatorFunctionNet, 537
RoboError_GilsonWrondID	Send
CRoboDeviceNet, 403	CSerialPortNet, 437
RoboError_McsBus_UnknownCommand	SendBuffered
CRoboDeviceNet, 403	CGilsonDeviceNet, 101
RoboError_NoEndSwitch	SendChannelData
CRoboDeviceNet, 403	CDigOutStimulatorFunctionNet, 62
RoboError_NoMoreData	CStg200xDownloadBasicNet, 474
CRoboDeviceNet, 403	SendCommand
RoboError_NoReference	CLIH3DeviceNet, 126
CRoboDeviceNet, 403	SendImmediate
RoboError NoSpeedOrAcceleration	CGilsonDeviceNet, 101
CRoboDeviceNet, 403	SendImmediateGetResponse
RoboError_OverPressure	CGilsonDeviceNet, 101
CRoboDeviceNet, 404	SendMultiplexedData
RoboError ParameterNotAllowed	CStimulusFunctionNet, 492
CRoboDeviceNet, 404	SendPreparedData
RoboError_PeristalticTimeout	CStimulusFunctionNet, 492
CRoboDeviceNet, 404	CW2100_StimulatorFunctionNet, 537
RoboError_Phase0OutOfRange	SendSegmentDefine
CRoboDeviceNet, 404	CStg200xDownloadNet, 482
RoboError_PollLoopCanceled	SendSegmentSelect
CRoboDeviceNet, 404	CStg200xDownloadNet, 482
RoboError_PollLoopCanceledAndStopMovement	SendSegmentStart
CRoboDeviceNet, 404	-
	CStg200xDownloadNet, 483 SendStart
RoboError_Pressure	
CRoboDeviceNet, 404	CStg200xBasicNet, 454
RoboError_RangeExceeded	CStimulusFunctionNet, 492
CRoboDeviceNet, 404	CW2100_StimulatorFunctionNet, 537
RoboError_StateChangeNotPossible	SendStartDacq
CRoboDeviceNet, 404	CMcsUsbDacqNet, 202
RoboError_Timeout	SendStartStgAndDacq
CRoboDeviceNet, 404	CMcsUsbDacqNet, 202
RoboError_UnknownCommand	SendStop
CRoboDeviceNet, 405	CStg200xBasicNet, 454
RoboMainLowLevelCommand	CStimulusFunctionNet, 492
CRoboDeviceNet, 405	CW2100_StimulatorFunctionNet, 538
RoboMainStatorLowLevelCommand	SendStopDacq
CRoboStatorDeviceNet, 417	CMcsUsbDacqNet, 203
RoboStatusEvent	SendStopStgAndDacq
CRoboDeviceNet, 405	CMcsUsbDacqNet, 203
RoboStatusEventDelegate	SendStopStgAndDacqWithOptions
Mcs::Usb, 29	CMcsUsbDacqNet, 203
RunTable	SendSyncData
CRoboDacqNet, 386	CStg200xDownloadBasicNet, 474
Eth Au tha Ala	Sensor
s_FilterAttributesNet	CFYIDeviceNet, 81
CFilterCoefficientsNet::s_FilterAttributesNet, 615	CMeasureTableDeviceNet, 292
SampleRate	CPatchServerDeviceNet, 326
CCreateFilterNet, 50	SerialNumber
Samplerate	CMcsUsbListEntryNet, 236
CMcsUsbDacqNet, 219	CMcsUsbNet, 261
Scale	SerialPort
CCreateFilterNet, 50	

CHLADeviceNet, 112	SetAnalogOutADCRange
Set4ADCCatchampAverageShift	CSCUFunctionNet, 433
CMcsBus_SensorNet, 158	SetAnalogOutChannel
Set4ADCMode	CW2100_FunctionNet, 530
CMcsBus SensorNet, 158	SetAnalogOutChannels
Set4DAC	CSCUFunctionNet, 433
CMcsBus_SensorNet, 158	SetAnalogOutDACRange
Set_Values	CSCUFunctionNet, 435
CNF GenDeviceNet, 317	SetAnalogOutFilter
	_
CPathIdentDeviceNet, 327	CW2100_FunctionNet, 530
SetAbsMaxCurrentInMicroAmp	SetAnalogThresholdHigh
CMultiwellOptoStimFunctionNet, 315	CWarnerValveControllerDeviceNet, 572
SetAccelerationI	SetAnalogThresholdLow
CRoboStatorDeviceNet, 414	CWarnerValveControllerDeviceNet, 572
SetAccelerationNativel	SetAnalogVoltageRange
CRoboStatorDeviceNet, 414	CPPCFunctionNet, 352
SetAccelerationNativeXY	SetAnalogVoltages
CRoboStatorDeviceNet, 414	CPPS_FunctionNet, 357
SetAccelerationNativeZ	SetAttenuation
CRoboStatorDeviceNet, 414	CChannelTestDeviceNet, 33
SetAccelerationXY	SetAudioChannels
CRoboStatorDeviceNet, 415	CMeaAudioFunctionNet, 265, 266
SetAccelerationZ	CW2100_FunctionNet, 530
CRoboStatorDeviceNet, 415	SetAudioOutDacParameter
SetAccelGyroDesiredRate	CLIH3DeviceNet, 127
CW2100_FunctionNet, 529	SetAutocalibrationDisabled
SetAccelGyroEnabled	CStg200xBasicNet, 454
CW2100_FunctionNet, 529	SetAxisConfig
	CRoboDeviceNet::RoboMainLowLevelCommands,
SetAccelRange	
CW2100_FunctionNet, 529	612
SetActiveRunningTableNumber	SetAxisLED
CWarnerValveControllerDeviceNet, 571	CRoboocyte2DeviceNet, 411
SetADC	SetAxisParametersEeprom
CWarnerValveControllerDeviceTesterFunctionNet,	CMcsBus_AxisParametersNet, 130
582	SetBandwidthByIndex
SetAdcChannels	CIntanMea_FunctionNet, 116
CSafeISDeviceNet, 418	SetBaseFrequency
SetADCInputOffset	CRFFunctionNet, 373
CCMOSMea_FunctionNet, 43	SetBaseSamplerate
SetAdcOffset	CCMOSMeaDeviceNet, 47
CLIH3DeviceNet, 126	SetBath
COctoPotDeviceNet, 320	CCMOSMea_FunctionNet, 43
SetAdcOffsetPermanent	SetBathclamp
CLIH3DeviceNet, 127	COctoPotDeviceNet, 320
SetAdcSamplePos	SetBathMode
CSafeISDeviceNet, 418	CCMOSMea_FunctionNet, 43
SetAirpressureLimit	SetBlankingEnable
CRoboDeviceNet, 398	CStg200xBasicNet, 455, 456
SetAirValve	SetBoostAlwaysOnMode
	•
CRoboDeviceNet, 399	CW2100_StimulatorFunctionNet, 538
SetAllDigout	SetBuffer
CRoboDacqNet, 386	CGenericDevelopDeviceNet, 96
SetAmplificationSwitch	0 ID (())
00 - 0 - 0 - 1 - 11 - 05	SetBufferIndex
COctoPotDeviceNet, 320	CTEERFunctionNet, 519
SetAmplitude	CTEERFunctionNet, 519 SetBusAddress
SetAmplitude CChannelTestDeviceNet, 33	CTEERFunctionNet, 519 SetBusAddress CMcsBusNet, 170
SetAmplitude	CTEERFunctionNet, 519 SetBusAddress

SetByteBuffer	CRoboDacqNet, 387
CGenericDevelopDeviceNet, 96	SetConfigurationBitRelais
SetCalibration	CRoboDacqNet, 387
CTcxDeviceNet, 508	SetConfigurationBitRV_Gen
SetCardinalDacqSamplerate	CRoboDacqNet, 388
CInterfaceboardFunctionNet, 119	SetConfigurationBitStream
SetCardinalStgOutputrate	CRoboDacqNet, 388
CInterfaceboardFunctionNet, 120	SetConfigurationBitSupply
SetChannel	CRoboDacqNet, 388
CSw2to64DeviceNet, 497	SetControllerParams
SetChannelmap	CTEERFunctionNet, 519
CWClassicFunctionNet, 586	SetCrossTalkOffset
SetChannels	CRoboDacqNet, 388
CSw2to64DeviceNet, 497	SetCrossTalkOptimum
SetChannelSwitch	CRoboDacqNet, 388
COctoPotDeviceNet, 320	SetCurrentAirvalveLimit
	CRoboDeviceNet, 399
SetChargingMode	•
CMultiBatteryChargerDeviceNet, 302	SetCurrentAndAir
SetChargingPCoefficient	CRoboDeviceNet, 399
CMultiBatteryChargerDeviceNet, 302	SetCurrentAndAirXY
SetCheckVoltage	CRoboStatorDeviceNet, 415
COkuvisionStimulatorDeviceNet, 324	SetCurrentEditTableNumber
SetClampMode	CWarnerValveControllerDeviceNet, 572
CTEERFunctionNet, 519	SetCurrentEnable
CWarnerUssingFunctionNet, 553	CTEERFunctionNet, 519
SetColorRgb	SetCurrentFactor
CMultiwellOptoStimFunctionNet, 315	COkuvisionStimulatorDeviceNet, 324
SetColorStr	SetCurrentMode
CMultiwellOptoStimFunctionNet, 316	CStg200xBasicNet, 456
SetCommand	SetCycles
	SetCycles CMeaCleanDeviceNet, 269
SetCommand CMcsBusNet, 170, 171 CPedoterDeviceNet, 328	-
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328	CMeaCleanDeviceNet, 269
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode CSafeISDeviceNet, 419
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode CSafeISDeviceNet, 419
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode CSafeISDeviceNet, 419 SetDacPulseform
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode CSafeISDeviceNet, 419 SetDacPulseform CSafeISDeviceNet, 419
CMcsBusNet, 170, 171 CPedoterDeviceNet, 328 CRoboDacqNet, 386 SetConfiguration CMcsUsbNet, 254 SetConfigurationBit CRoboDacqNet, 386 SetConfigurationBitAxc CRoboDacqNet, 386 SetConfigurationBitBlu_Led CRoboDacqNet, 386 SetConfigurationBitBlu_LedToggleFast CRoboDacqNet, 387 SetConfigurationBitBlu_LedToggleSlow CRoboDacqNet, 387 SetConfigurationBitCC_Gen CRoboDacqNet, 387 SetConfigurationBitCV_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRC_Gen CRoboDacqNet, 387 SetConfigurationBitRed_Led CRoboDacqNet, 387 SetConfigurationBitRed_LedSaturation CRoboDacqNet, 387 SetConfigurationBitRed_LedSaturation CRoboDacqNet, 387 SetConfigurationBitRed_LedToggleFast	CMeaCleanDeviceNet, 269 CMeaCoatDeviceNet, 274 SetD CTcxDeviceNet, 509 SetDacAmplificationFactor CStg200xBasicNet, 456 SetDacAutoControl COctoPotDeviceNet, 320 SetDacIdleValue CLIH3DeviceNet, 127 SetDacMode CSafeISDeviceNet, 418 SetDACOffset CGrapheneFunctionNet, 106 COkuvisionStimulatorDeviceNet, 324 SetDacOffset CDacCalibrationFunctionNet, 51 CLIH3DeviceNet, 127 COctoPotDeviceNet, 320 SetDacOffsetPermanent CLIH3DeviceNet, 128 SetDacPeriode CSafeISDeviceNet, 419 SetDacPulseform CSafeISDeviceNet, 419 SetDacQLegacyMode

CW2100_FunctionNet, 530	SetDownsampleFactor
SetDACs	CRoboDacqNet, 388
CMcsBus_SensorNet, 158	SetDSPHighPassByIndex
SetDacUseIdleValue	CIntanMea_FunctionNet, 116
CLIH3DeviceNet, 128	SetDuration
SetDacValue	CMeaCoatDeviceNet, 274
COctoPotDeviceNet, 321	SetEEpromPage
SetDataMode	CLIH3DeviceNet, 128
CMcsUsbDacqNet, 204	SetElectrodeDacMux
SetDefault	CStg200xBasicNet, 457-459
CWarnerValveControllerDeviceNet, 572	SetElectrodeEnable
SetDestinationSerialNumber	CStg200xBasicNet, 459-461
CMcsUsbFactoryNet, 227	SetElectrodeMode
SetDetectionThreshold	CStg200xBasicNet, 462, 463
CMcsBus_SensorNet, 159	SetEnableAmplifierProtectionSwitch
SetDevice	CStg200xBasicNet, 463, 464
CTcxDeviceNet, 509	SetEnableHeaterLimit
SetDeviceList	CTcxDeviceNet, 509
CPositionImpDeviceNet, 343	SetEnablePulse
SetDeviceType	CWarnerUssingFunctionNet, 553
CTcxDeviceNet, 509	SetEnableThermocouple
SetDevname	CTcxDeviceNet, 509
CTcxDeviceNet, 509	SetExternalElectrodeEnable
SetDiagnosticMode	CStg200xBasicNet, 465
CIntanMea_FunctionNet, 116	SetExternalLED
SetDigitalData CMacDigitalDataFunctionNet 285	CTEERFunctionNet, 520
CMeaDigitalDataFunctionNet, 285	SetFAAmplification
SetDigitalOut	CStg200xBasicNet, 466
CMeaDeviceNet, 280	SetFilter OP 1 P 2 P 2 P 2 P 2 P 2 P 2 P 2 P 2 P 2
SetDigitalOutPortInvert	CRoboDacqNet, 388
CWarnerValveControllerDeviceNet, 572	SetFilterCoeffs
SetDigitalOutPortValve	CRoboDacqNet, 389
CWarnerValveControllerDeviceNet, 573	SetFilterParameter
SetDigitalPortDirection	CFilterConfigurationNet, 69
CWarnerValveControllerDeviceNet, 573	CFilterConfigurationRegisterNet, 71
SetDigitalSource	SetFilterParameterPermanent
CMcsUsbDacqNet, 204–206	CFilterConfigurationNet, 69
SetDigitalStimulatorTrigger	CFilterConfigurationRegisterNet, 71
CW2100_StimulatorFunctionNet, 538	SetFilterParametersHeadstage
SetDigitalStimulatorTriggerSlope	CWClassicFunctionNet, 586
CW2100_StimulatorFunctionNet, 538	SetFinalDischargeVoltage
SetDigout	CMultiBatteryChargerDeviceNet, 303
CFluidControlDeviceNet, 78	SetFrequency
CRoboDacqNet, 388	CChannelTestDeviceNet, 33
SetDigoutMode	CRadioControledDevicesNet, 366
CStg200xBasicNet, 457	SetGain
SetDigOutState	CPgaDeviceNet, 332
CLIH3DeviceNet, 128	SetGate
SetDigoutValue	CCMOSMea_FunctionNet, 43
CStg200xBasicNet, 457	SetGateFloating
SetDIO	CCMOSMea_FunctionNet, 44
CMcsBus_FYIExtensionNet, 132	SetGateToVOP
SetDischargeCurrentSetPoint	CCMOSMea_FunctionNet, 44
CMultiBatteryChargerDeviceNet, 303	SetGlobalRepeat
SetDisplayMode	CDigOutStimulatorFunctionNet, 62
CWarnerValveControllerDeviceNet, 573	SetGyroRange
SetDisplayText	CW2100_FunctionNet, 530
CRoboDacqNet, 388	SetHasChecksum

CWClassicFunctionNet, 586	SetLatency
SetHeadstage	CMcsBus_SensorNet, 159
CStg200xBasicNet, 466	SetLayoutConfiguration
SetHeadstageOnOff	CMEA2100x256FunctionNet, 262
CW2100_FunctionNet, 530	SetLED
CWClassicFunctionNet, 586	CRetinaLedDeviceNet, 368
SetHeadstageSamplingActive	SetLEDSwitch
CW2100_FunctionNet, 530	CMcsBus_ExtensionNet, 131
SetHeadstageToSleep	SetLength
CW2100_FunctionNet, 530	CRobo_FYIProgram_FunctionNet, 375
SetHeaterLimit	SetLiquidResistance
CTcxDeviceNet, 509	CTEERFunctionNet, 520
SetHighCurrentMode	CWarnerUssingFunctionNet, 554
CWarnerUssingFunctionNet, 553	SetListmodeIndexRange
SetHighpassFilterEnable	CStg200xBasicNet, 466
CFilterConfigurationNet, 70	SetListmodeTriggerSource
SetHWConfig	CStg200xBasicNet, 466
CRoboDeviceNet::RoboMainLowLevelCommands,	SetLowCurrentMode
612	CWarnerUssingFunctionNet, 554
SetHWRevision	SetLumi
CRoboDeviceNet::RoboMainLowLevelCommands,	CRetinaLedDeviceNet, 368
612	SetMaxCurrent
SetHWRevisionEeprom	CMeaCoatDeviceNet, 275
CMcsBusNet, 171	SetMaxDurationHighCurrentInMicroSec
SetHWSelectedChannels	CMultiwellOptoStimFunctionNet, 316
CWClassicFunctionNet, 586	SetMaxDutyCycleHighCurrent
Setl	CMultiwellOptoStimFunctionNet, 316
CTcxDeviceNet, 509	SetMaxHeaterPowerMultiwell
	CTcxDeviceNet, 510
SetIClamp	
CRoboDacqNet, 389	SetMaxNoPressure
SetlCoeff CDaha EVITama FunctionNet 277	CRoboDeviceNet::RoboMainLowLevelCommands,
CRobo_FYITemp_FunctionNet, 377	612
SetICOffset	SetMaxNoPressureWaitTime
CRoboDacqNet, 389	CRoboDeviceNet::RoboMainLowLevelCommands,
SetIdleModeOffset	612
CWarnerUssingFunctionNet, 554	SetMaxP
SetIGain	CTcxDeviceNet, 510
CRoboDacqNet, 389	SetMaxPower
SetImpedanceTestFrequency	COkuvisionStimulatorDeviceNet, 325
CMealmpedanceDeviceNet, 291	CRobo_FYITemp_FunctionNet, 377
SetImpId	SetMaxPressureWaitTime
CPositionImpDeviceNet, 343	${\tt CRoboDeviceNet::} RoboMainLowLevel Commands,$
SetImplantCurrentSetpoint	612
CPositionIIDeviceNet, 338	SetMaxVoltage
SetInMovement	CMeaCleanDeviceNet, 269
CRoboDeviceNet, 399	COkuvisionStimulatorDeviceNet, 325
SetIntanRegister	SetMCAcceleration
CIntanMea_FunctionNet, 116	CMcsBus_MotorControlNet, 143
SetIntBuffer	SetMCAccelerationEeprom
CGenericDevelopDeviceNet, 96	CMcsBus_MotorControlNet, 143
SetIO	SetMCAccelerationShortCommand
CWarnerValveControllerDeviceTesterFunctionNet,	CMcsBus_MotorControlNet, 143
582	SetMCAxisRevisionEeprom
SetIODirection	CMcsBus_MotorControlNet, 143
CWarnerValveControllerDeviceTesterFunctionNet,	SetMCBreakCurrent
582	CMcsBus_MotorControlNet, 143
SetloVoltage	SetMCBreakCurrentEeprom
CInterfaceboard2FunctionNet, 118	CMcsBus MotorControlNet. 144

SetMCConfig	SetMCSpeedEeprom
CMcsBus_MotorControlNet, 144	CMcsBus_MotorControlNet, 149
SetMCConfigEeprom	SetMCSpeedShortCommand
CMcsBus_MotorControlNet, 144	CMcsBus_MotorControlNet, 149
SetMCCurrent	SetMCSpeedUnitEeprom
CMcsBus_MotorControlNet, 144	CMcsBus_MotorControlNet, 149
SetMCCurrentEeprom	SetMCStandbyCurrent
CMcsBus_MotorControlNet, 144	CMcsBus_MotorControlNet, 149
SetMCCurrentMode	SetMCStandbyCurrentEeprom
CMcsBus MotorControlNet, 144	CMcsBus_MotorControlNet, 149
SetMCCurrentModeEeprom	SetMCStandbyTime
CMcsBus_MotorControlNet, 145	CMcsBus_MotorControlNet, 149
SetMCCurrentModeShortCommand	SetMCStandbyTimeEeprom
CMcsBus_MotorControlNet, 145	CMcsBus_MotorControlNet, 150
SetMCCurrentPosition	SetMeasurementMode
CMcsBus_MotorControlNet, 145	CStg200xBasicNet, 466
SetMCCurrentShortCommand	SetMinimalThreshold
CMcsBus_MotorControlNet, 145	CMcsBus SensorNet, 159
SetMCMaxAcceleration	SetMinNoPressureWaitTime
CMcsBus_MotorControlNet, 145	CRoboDeviceNet::RoboMainLowLevelCommands
SetMCMaxAccelerationEeprom	612
CMcsBus_MotorControlNet, 145	SetMinPressure
SetMCMaxCurrent	CRoboDeviceNet, 399
CMcsBus_MotorControlNet, 146	CRoboDeviceNet::RoboMainLowLevelCommands
SetMCMaxCurrentEeprom	612
CMcsBus_MotorControlNet, 146	SetMinPressureWaitTime
SetMCMaxSpeed	CRoboDeviceNet::RoboMainLowLevelCommands
CMcsBus_MotorControlNet, 146	612
SetMCMaxSpeedEeprom	SetMinVoltage
CMcsBus_MotorControlNet, 146	CMeaCleanDeviceNet, 270
SetMCMaxTravel	SetModeSelect
CMcsBus_MotorControlNet, 146	CPulseGeneratorFunctionNet, 364
SetMCMaxTravelEeprom	SetMovePump
CMcsBus_MotorControlNet, 146	CMcsBus SensorNet, 159
SetMCMaxTravelShortCommand	SetMultiHeadstageMode
	-
CMcsBus_MotorControlNet, 147 SetMCNewPosition	CW2100_FunctionNet, 531 SetNanoVoltsPerKelvin
CMcsBus_MotorControlNet, 147	CMcsBus_TempSensorNet, 162
SetMCOutputOnOff CMapRya MaterControlNet 147	SetNeurochipMemoryData
CMcsBus_MotorControlNet, 147	CCMOSMea_FunctionNet, 44
SetMCReference CMacRus MaterControlNet 147	SetNoFilterCoeffs CRaha Reserved. 200
CMcsBus_MotorControlNet, 147	CRoboDacqNet, 389
SetMCReferenceCurrent CMasRus MaterControlNet 147	SetNumberOfAnalogChannels
CMcsBus_MotorControlNet, 147	CMeaDeviceNet, 280
SetMCReferenceCurrentEeprom	SetNumberOfChannels
CMcsBus_MotorControlNet, 147	CMeaDeviceNet, 281, 282
SetMCRegulatorGain	COctoPotDeviceNet, 321
CMcsBus_MotorControlNet, 148	SetOffsetCurrent
SetMCRegulatorGainEeprom	CMeaCoatDeviceNet, 275
CMcsBus_MotorControlNet, 148	SetOnOff
SetMCRotation	CTcxDeviceNet, 510
CMcsBus_MotorControlNet, 148	SetOutputMap
SetMCScalingFactor	CStg200xDownloadNet, 483
CMcsBus_MotorControlNet, 148	SetOutputRate
SetMCScalingFactorEeprom	COctoPotDeviceNet, 321
CMcsBus_MotorControlNet, 148	CStg200xBasicNet, 467
SetMCSpeed	SetP
CMcsBus_MotorControlNet, 148	CTcxDeviceNet, 510

SetParameter	SetPumpManualOnOff
CRoboDeviceNet::RoboMainLowLevelCommands,	CPPS_FunctionNet, 358
613	SetPumpMaxSpeed
SetPattern	CPPS FunctionNet, 358
CMeaSwitchDeviceNet, 294	SetPumpModeType
SetPatternBool	CPPCFunctionNet, 353
CMeaSwitchDeviceNet, 294	CPPS_FunctionNet, 358
SetPatternListEntry	SetPumpSpeed
COctoPotDeviceNet, 321	CRoboFluidDeviceNet, 407
SetPauseDuration	SetPumpSpeedRatio
	• •
CMeaCoatDeviceNet, 275	CPPS_FunctionNet, 359
SetPCoeff	SetPumpSpeedUnit
CRobo_FYITemp_FunctionNet, 377	CPPCFunctionNet, 353
SetPeriod	CPPS_FunctionNet, 359
CPulseGeneratorFunctionNet, 364	SetPWM
SetPeriod_us	CFluidControlDeviceNet, 78
CTEERFunctionNet, 520	SetRampParameter
SetPermanentCurrentInMicroAmp	COctoPotDeviceNet, 321
CMultiwellOptoStimFunctionNet, 316	SetRatedCapacity
SetPersistency	CMultiBatteryChargerDeviceNet, 303
CRetinaLedDeviceNet, 368	SetRatedCapacityVolatile
SetPGain	CMultiBatteryChargerDeviceNet, 303
CRoboDacqNet, 389	SetRecordingNumber
SetPidParameter	CRoboDacqNet, 389
COctoPotDeviceNet, 321	SetReferenceElectrodeMode
SetPiezoState	CSCUFunctionNet, 435
CMcsBus_SensorNet, 159	SetReferenceElectrodeSwitchState
SetPlateMux	CSCUFunctionNet, 435
CMultiwellDeviceNet, 310, 311	SetRegionOfInterests
SetPlateType	CCMOSMeaDeviceNet, 47
CMultiwellDeviceNet, 311	SetRegulationTimeouts
SetPoti	CMcsBus_SensorNet, 159
CMcsUsbDacqNet, 206	SetRegulatorFactor
SetPowerMuxPlate	CMcsBus_SensorNet, 160
CMultiwellDeviceNet, 311	SetRegulatorOnOff
SetPowerStrength	CMcsBus_SensorNet, 160
CPositionIIDeviceNet, 338	CRobo_FYITemp_FunctionNet, 377
SetPressureOffset	SetRepeat
CMcsBus_SensorNet, 159	CRetinaLedDeviceNet, 369
CPPCFunctionNet, 352	SetRepeats
SetPressureRange	CProgramPressureCurveNet, 361
CPPCFunctionNet, 352	SetResetFilter
SetPulse	CWClassicFunctionNet, 586
CWarnerUssingFunctionNet, 554	
•	SetRFFrequency
SetPulseform	CPositionImpDeviceNet, 344
COkuvisionStimulatorDeviceNet, 325	SetRFFrequencyHeadstage
SetPulseLength	CWClassicFunctionNet, 586
CPulseGeneratorFunctionNet, 364	SetRFFrequencyReceiver
SetPumpCouple	CWClassicFunctionNet, 587
CPPS_FunctionNet, 358	SetRFFrequencyReceiverEeprom
SetPumpEnableSpeedRatio	CWClassicFunctionNet, 587
CPPS_FunctionNet, 358	SetRFLostBehaviour
SetPumpFastOnOff	CWClassicFunctionNet, 587
CPPS_FunctionNet, 358	SetRFPower
SetPumpFastSpeed	CWClassicFunctionNet, 587
CPPS_FunctionNet, 358	SetRotatePump
SetPumpFunctionSpeeds	CMcsBus_SensorNet, 160
CPPS_FunctionNet, 358	SetRTC
Of 1 O_1 unotioningt, 000	OGII II O

COkuvisionStimulatorDeviceNet, 325	CMeaCleanDeviceNet, 270
CPositionIIDeviceNet, 338	CMeaCoatDeviceNet, 275
SetSampleInterval	SetSoftwareKey
CLIH3DeviceNet, 129	CMcsUsbNet, 254
SetSamplePeriode	SetSollPressure
CMcsBus_SensorNet, 160	CMcsBus_SensorNet, 160
SetSamplerate	SetSollTemp
CMcsUsbDacqNet, 206	CRobo_FYITemp_FunctionNet, 377
SetScreen	SetSourceBulk
CRoboDacqNet, 389	CCMOSMea_FunctionNet, 44
SetSearchReferenceFastAccel	SetSourceDrain
CRoboDeviceNet::RoboMainLowLevelCommands,	CCMOSMea_FunctionNet, 44
613	SetSourceGate
SetSearchReferenceFastSpeed	CCMOSMea_FunctionNet, 44
CRoboDeviceNet::RoboMainLowLevelCommands,	SetSpeedI
613	CRoboStatorDeviceNet, 415
SetSearchReferenceFineAccel	SetSpeedNativel
CRoboDeviceNet::RoboMainLowLevelCommands,	CRoboStatorDeviceNet, 415
613	SetSpeedNativeXY
SetSearchReferenceFineSpeed	CRoboStatorDeviceNet, 415
CRoboDeviceNet::RoboMainLowLevelCommands,	SetSpeedNativeZ
613	CRoboStatorDeviceNet, 415
SetSearchReferenceMethod	SetSpeedXY
CRoboDeviceNet::RoboMainLowLevelCommands,	CRoboStatorDeviceNet, 415
613	SetSpeedZ
SetSearchReferenceMoveOut	CRoboStatorDeviceNet, 416
CRoboDeviceNet::RoboMainLowLevelCommands,	SetStartTriggerSlope
614	
	CDigOutStimulatorFunctionNet, 62
SetSearchReferenceOffsetPos	SetStateDebugData
CRoboDeviceNet::RoboMainLowLevelCommands,	CPositionIIDeviceNet, 340
614	SetStateEventData
SetSelectedChannels	CPositionIIDeviceNet, 340
CMcsUsbDacqNet, 207–209	SetStgProgramInfo
CW2100_FunctionNet, 531	CStg200xBasicNet, 467
SetSelectedChannelsQueue	SetStimulusSites
CMcsUsbDacqNet, 209-211	CCMOSMea_FunctionNet, 44
SetSelectedData	SetStopTriggerSlope
CMcsUsbDacqNet, 211-213	CDigOutStimulatorFunctionNet, 63
SetSelectedHeadstage	SetStringFormat
CWClassicFunctionNet, 587	CMcsUsbListEntryNet, 235
SetSensorType	CMcsUsbListNet, 239
CTcxDeviceNet, 510	SetSubChannel
SetSerialNumberHeadstage	CMcsBus_MotorControlNet, 150
CWClassicFunctionNet, 587	SetSwitches
	CSafeISDeviceNet, 419
SetSetpoint Standard 510	•
CTcxDeviceNet, 510	SetSyncoutMap
SetShortBuffer	CStg200xBasicNet, 467
CGenericDevelopDeviceNet, 97	SetTableName
SetSimulation	CWarnerValveControllerDeviceNet, 573
CRoboDacqNet, 389	SetTablepointer
SetSineParameter	CRetinaLedDeviceNet, 369
COctoPotDeviceNet, 321	SetTableStep
SetSingleHeater	CWarnerValveControllerDeviceNet, 574
CMcsBus_FYIExtensionNet, 133	SetTableStepAll
SetSingleValve	CWarnerValveControllerDeviceNet, 574
CFluidControlDeviceNet, 78	SetTestMode
CRoboFluidDeviceNet, 408	CRFFunctionNet, 373
SetSlope	SetThermocoupleNanovoltPerKelvin
In a	

CFluidControlDeviceNet, 80	SetValveDigitalInInvert
CTcxDeviceNet, 511	CWarnerValveControllerDeviceNet, 574
SetThermoOffset	SetValveDigitaIInPort
CMcsBus_TempSensorNet, 162	CWarnerValveControllerDeviceNet, 575
SetTransformer	SetValveLedOn
CMeFunctionNet, 296	CWarnerValveControllerDeviceNet, 575
SetTrigger	SetValveManualGroup
CRetinaLedDeviceNet, 369	CWarnerValveControllerDeviceNet, 575
CWarnerValveControllerDeviceTesterFunctionNet,	SetValveManualState
582	CWarnerValveControllerDeviceNet, 575
SetTriggerMaskValue	SetValveMode
CMeaDeviceNet, 282	CWarnerValveControllerDeviceNet, 576
SetTriggerPeriod	SetValves
CMeaDeviceNet, 283	CMcsBus_FYIExtensionNet, 133
SetTriggerSource	SetValvesActiveMap
CStg200xBasicNet, 467, 468	CWarnerValveControllerDeviceNet, 576
SetTriggerSyncDirection	SetValvesManualStateMap
CWarnerValveControllerDeviceTesterFunctionNet,	CWarnerValveControllerDeviceNet, 576
583	SetValveTableEntry
SetUByteBuffer	CWarnerValveControllerDeviceNet, 576
CGenericDevelopDeviceNet, 97	SetVds
SetUClamp	CGrapheneFunctionNet, 107
CRoboDacqNet, 390	SetVdVs
SetUCOffset	CGrapheneFunctionNet, 107
CRoboDacqNet, 390	SetVdVsDAC
SetUIntBuffer	CGrapheneFunctionNet, 108
CGenericDevelopDeviceNet, 98	SetVelocityI
SetupGroupDacqQueue	CRoboStatorDeviceNet, 416
CMcsUsbDacqNet, 213	SetVelocityXY
SetupRetriggerMode	CRoboStatorDeviceNet, 416
CStg200xDownloadBasicNet, 475	SetVelocityZ
SetupTrigger	CRoboStatorDeviceNet, 416
CStg200xDownloadBasicNet, 475	SetVgs
CStimulusFunctionNet, 493	CGrapheneFunctionNet, 108
SetupTriggerSingle	SetVMMaxNegativeCurrent
CStg200xDownloadBasicNet, 476	CMcsBus_VoltageModeNet, 165
CStimulusFunctionNet, 493	SetVMMaxNegativeCurrentEeprom
SetUseBubble	CMcsBus_VoltageModeNet, 166
CPPS_FunctionNet, 359	SetVMMaxNegativeVoltage
SetUserParameter	CMcsBus VoltageModeNet, 166
CRoboDeviceNet::RoboMainLowLevelCommands,	SetVMMaxNegativeVoltageEeprom
614	CMcsBus_VoltageModeNet, 166
SetUShortBuffer	SetVMMaxPositiveCurrent
CGenericDevelopDeviceNet, 99	CMcsBus_VoltageModeNet, 166
SetUVOffset	SetVMMaxPositiveCurrentEeprom
CRoboDacqNet, 390	CMcsBus_VoltageModeNet, 166
SetValue	SetVMMaxPositiveVoltage
CGenericDevelopDeviceNet, 99	CMcsBus VoltageModeNet, 166
SetValve	SetVMMaxPositiveVoltageEeprom
CFluidControlDeviceNet, 80	CMcsBus_VoltageModeNet, 167
CRoboFluidDeviceNet, 408	SetVMOutputOnOff
SetValve1	CMcsBus_VoltageModeNet, 167
CRobo_FYIProgram_FunctionNet, 375	SetVMVoltage
SetValve2	CMcsBus_VoltageModeNet, 167
CRobo_FYIProgram_FunctionNet, 375	SetVoltage12VLimit
SetValveActive	CRoboDeviceNet, 399
CPPCFunctionNet, 353	SetVoltage5VLimit
CWarnerValveControllerDeviceNet. 574	CRoboDeviceNet, 399

SetVoltageAirvalveLimit	CMcsUsbDacqNet, 216, 217
CRoboDeviceNet, 400	StartMCMovement
SetVoltageClampControllerParam_D	CMcsBus_MotorControlNet, 150
CWarnerUssingFunctionNet, 555	StartMeasurement
SetVoltageClampControllerParam_I	CMealmpedanceDeviceNet, 291
CWarnerUssingFunctionNet, 555	StartPoll
SetVoltageClampControllerParam_P	CStimulusFunctionNet, 495
CWarnerUssingFunctionNet, 555	CW2100 StimulatorFunctionNet, 538
SetVoltageMode	StartSampling
CStg200xBasicNet, 468	CTEERFunctionNet, 521
SetVoltageRange	StartSync
CGrapheneFunctionNet, 109	CMcsBus_SensorNet, 160
SetVoltageRangeByIndex	State
CMcsUsbDacqNet, 213	HeadStageIDTypeState, 604
SetVoltageRangeInMicroVolt	Status
CMcsUsbDacqNet, 214	CUsbExceptionNet, 523
SetVoltageResolution	Status AlreadyConfigured
CGrapheneFunctionNet, 109	CMcsUsbNet, 256
SetVoltageRs485ALimit	Status BadStartFrame
CRoboDeviceNet, 400	CMcsUsbNet, 257
SetVoltageRs485BLimit	Status Btstuff
CRoboDeviceNet, 400	CMcsUsbNet, 257
SetVoltageValvesLimit	Status BufferOverrun
CRoboDeviceNet, 400	CMcsUsbNet, 257
SetWaveform	Status BufferUnderrun
CChannelTestDeviceNet, 34	CMcsUsbNet, 257
CTEERFunctionNet, 520	Status_Canceled
SetWaveLengthInNanometer	CMcsUsbNet, 257
CMultiwellOptoStimFunctionNet, 317	Status_Canceling
SetWorkingFrequency	CMcsUsbNet, 257
CRFFunctionNet, 373	Status_ConnectedPipes
SetWPADebugMode	CMcsUsbNet, 257
CWClassicFunctionNet, 587	Status_ControlNotOwned
SetWPAType	CMcsUsbNet, 257
CWClassicFunctionNet, 587	Status_Crc
SetXGain	CMcsUsbNet, 257
	,
CRoboDacqNet, 390	Status_DataOverrun
Sideband	CMcsUsbNet, 257
CStimulusFunctionNet::SidebandData, 617	Status_DataToggleMismatch
SidebandData	CMcsUsbNet, 257
CStimulusFunctionNet::SidebandData, 617	Status_DataUnderrun
SineStart 200	CMcsUsbNet, 258
COctoPotDeviceNet, 322	Status_DeviceLocked
size	CMcsUsbNet, 258
DigitalSource< digitalsourceenum >, 591	Status_DeviceNotFound
SN	CMcsUsbNet, 258
HeadStageIDType, 602	Status_DeviceRemoved
Source	CMcsUsbNet, 258
DigitalSource< digitalsourceenum >, 591	Status_DevNotResponding
Start	CMcsUsbNet, 258
CMeaCleanDeviceNet, 270	Status_EndpointHalted
CMeaCoatDeviceNet, 275	CMcsUsbNet, 258
CRobo_FYIProgram_FunctionNet, 375	Status_ErrorBusy
StartDacq	CMcsUsbNet, 258
CMcsUsbDacqNet, 214, 215	Status_ErrorShortTransfer
StartInternalCalibration	CMcsUsbNet, 258
CTEERFunctionNet, 520	Status_Fifo
StartLoop	CMcsUsbNet, 258

Status_FrameControlOwned	StimulusDeviceDataAndUnrolledData
CMcsUsbNet, 258	CStimulus Function Net:: Stimulus Device Data And Unrolled Data,
Status_InternalHcError	619
CMcsUsbNet, 258	StimulusFunction
Status_InvalidParameter	CLIH3DeviceNet, 129
CMcsUsbNet, 259	StimulusParameters
Status_InvalidPipeHandle	HeadStageIDType, 602
CMcsUsbNet, 259	Stop
Status_InvalidUrbFunction	CMeaCleanDeviceNet, 270
CMcsUsbNet, 259	CMeaCoatDeviceNet, 276
Status_IoPending	StopDacq
CMcsUsbNet, 259	CMcsUsbDacqNet, 218
Status_loTimeout	StopLoop
CMcsUsbNet, 259	CMcsUsbDacqNet, 218
Status_IsochRequestFailed	StopMCMovement
CMcsUsbNet, 259	CMcsBus_MotorControlNet, 150
Status_LastUsbErrorMismatch	StopMovement
CMcsUsbNet, 259	CRoboDeviceNet, 400
Status_NoBandwidth	StopMovementI
CMcsUsbNet, 259	CRoboStatorDeviceNet, 416
Status_NoMemory	StopMovementXY
CMcsUsbNet, 259	CRoboStatorDeviceNet, 416
Status_NoSuchDevice	StopMovementZ
CMcsUsbNet, 259	CRoboStatorDeviceNet, 416
Status_NotAccessed	StopPlateClamp
CMcsUsbNet, 259	CMultiwellDeviceNet, 312
Status_NotSupported	StopPoll
CMcsUsbNet, 260	CStimulusFunctionNet, 495
Status_PidCheckFailure	CW2100_StimulatorFunctionNet, 538
CMcsUsbNet, 260	StopSampling
Status_PipeNotLinked	CTEERFunctionNet, 521
CMcsUsbNet, 260	StopTable
Status_RequestFailed	CRoboDacqNet, 390
CMcsUsbNet, 260	StoreValveTable
Status_RequestMutexFailed	CWarnerValveControllerDeviceNet, 577
CMcsUsbNet, 260	SwitchOnOff
Status_RequestMutexTimeout	CPositionIIDeviceNet, 340
CMcsUsbNet, 260	SYNC_BIT0
Status_Stall	CW2100_StimulatorFunctionNet, 539
CMcsUsbNet, 260	SYNC_BIT1
Status_Unconfigured	CW2100_StimulatorFunctionNet, 539
CMcsUsbNet, 260	T11 W %
Status_UnexpectedPid	Table_Wait
CMcsUsbNet, 260	CRoboDacqNet, 390
Stg200xPollStatusEvent	TableDefBegin
CStg200xDownloadNet, 483	CRoboDacqNet, 390
StgStatusNet, 618	TableDefEnd
FromIntPtr, 618	CRoboDacqNet, 390
FromPtr, 618	TableEntryChangedEvent
ListOfChangedTriggers, 618	CWarnerValveControllerDeviceNet, 580
TiggerStatus, 618	TactSwitchGetState
StillConnected	CMcsBus_SensorNet, 160
CRadioControledDevicesNet, 366	TactSwitchSetDisplay
Stimulator	CMcsBus_SensorNet, 161
CW2100_FunctionNet, 531	TEERFunctionNet
Stimulus	CTEERMachineDeviceNet, 522
CCMOSMeaDeviceNet, 47	ThrowCUsbExceptionNetOnError
CStg200xDownloadBasicNet, 477	CMcsUsbFunctionNet, 231
	CMcsUsbNet, 254

TiggerStatus	HeadStageIDType, 602
StgStatusNet, 618	W.P.I.
TimeResolutionInNanoSeconds	Valid
W2100_StimulusParametersNet, 621	HeadStageIDType, 602
ТоСрр	ValidKey
CFilterCoefficientsNet::s_FilterAttributesNet, 615	CMcsUsbNet, 255
ToString	Voltage
CFilterPropertyNet, 72	BatteryState, 30
CMcsUsbListEntryNet, 236	VoltageRangeDisplayStringMilliVolt
HeadStageIDType, 601	CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet
HeadstageIDTypeObject, 603	523
TriggerStatus	VoltageRangeInMicroVolt
CMcsUsbDeviceStatePushFunctionNet, 220	CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet
CMcsUsbDeviceStatePushNet, 221	523
TxnGetSerialNumber	W2100_StimulusParametersNet, 621
CMcsUsbNet, 254	VoltageResolutionInMicroVolt
TxnSetSerialNumber	W2100_StimulusParametersNet, 621
CMcsUsbNet, 254	VoltageString
TxnTestMemoryReadAndCheck	BatteryState, 30
CMcsUsbNet, 254	VOPSTimerSetResetTimes
TxnTestMemoryWrite	CCMOSMea_FunctionNet, 44, 45
CMcsUsbNet, 255	
Туре	W16lsW14
HeadStageIDType, 602	HeadStageIDType, 602
TypeValue	W2100_FunctionNet
HeadStageIDType, 602	CMeaDeviceNet, 284
7,,	W2100_StimulusParametersNet, 620
Unknown	CurrentRangeInNanoAmp, 621
HeadStageIDType, 601	CurrentResolutionInNanoAmp, 621
unknown	DACResolution, 621
Mcs::Usb, 27	TimeResolutionInNanoSeconds, 621
UnlockPlateClamp	VoltageRangeInMicroVolt, 621
CMultiwellDeviceNet, 312	VoltageResolutionInMicroVolt, 621
UnrolledAmplitude	WaitForAllChambers
CStimulusFunctionNet::StimulusDeviceDataAndUnro	blledD&WarnerUssingFunctionNet, 556
619	WaitForChamber
UnrolledDuration	CWarnerUssingFunctionNet, 556
CStimulusFunctionNet::StimulusDeviceDataAndUnro	bl le/diData UssingFunction
619	CWarnerUssingDeviceNet, 541
UnrolledSync	WClassicFunctionNet
CStimulusFunctionNet::StimulusDeviceDataAndUnro	olledD&MeaDeviceNet, 284
619	wIndex
UpdateChannelBlock	usbSetupPacket_t, 620
CCMOSMeaDeviceNet, 47	wLength
UpdateDisplay	usbSetupPacket_t, 620
CRoboDacqNet, 390	WPAError_ScanningIsPending
UpdateFirmware	CMcsUsbNet, 260
CMcsUsbFactoryNet, 227–229	Write
UpdateTransistorVoltages	CExternDTesterDeviceNet, 65
CCMOSMea_FunctionNet, 44	Write2
USB	CExternDTesterDeviceNet, 66
FirmwareDestinationNames, 600	WriteEepromRegisterPreconfig
usbSetupPacket_t, 620	CMcsUsbNet, 255
bmRequestType, 620	WritePipe
	CGenericDevelopDeviceNet, 100
bRequest, 620	·
wIndex, 620	WriteRegister CMod Johnst 255, 256
wLength, 620	CMcsUsbNet, 255, 256
wValue, 620	WriteRegister32 CMcsUshNet 256
LICAL INTRACINAMA	CAVICSUSDINEL ZOD

WriteRegisterArray
CMcsUsbNet, 256
WriteRegisterTimeSlot
CMcsUsbNet, 256
WriteRegisterValue
CMcsUsbNet, 256
wValue
usbSetupPacket_t, 620