

McsUsbNet.dll Version 5.1.30

Multi Channel Systems MCS GmbH
Aspenhaustrasse 21
72770 Reutlingen
Germany
Fon +49-71 21-90 92 5 - 0
Fax +49-71 21-90 92 5 -11
info@multichannelsystems.com
www.multichannelsystems.com

Generated by Doxygen 1.9.1

1	MCSUSDNet.dll for MCS USB devices	1
	1.1 Introduction	
	1.2 System requirements	
	1.3 Connecting to an MCS device	2
2	Poevice Classes	2
	2.1 The MCS FluidControl Device	2
	2.1.1 Introduction	2
	2.1.2 Access to the FluidControl device	2
	2.2 MCS-USB-Sw2to64 device	3
3	Function Classes	3
4	Data ACQuisition (DACQ) Devices	4
5	The MCS Robo Device	5
	5.1 Introduction	5
6	STG200x & STG400x STimulus Generator	5
	6.1 Introduction	5
	6.2 Download mode	6
	6.2.1 Memory Layout and Trigger Setup	6
	6.3 Streaming mode	8
	6.3.1 Memory Layout and Trigger Setup	9
7	Namespace Index	11
	7.1 Namespace List	11
8	Hierarchical Index	11
	8.1 Class Hierarchy	11
9	Class Index	16
	9.1 Class List	16
1	0 Namespace Documentation	22
	10.1 Mcs Namespace Reference	22
	10.2 Mcs::Usb Namespace Reference	23
	10.2.1 Enumeration Type Documentation	52
	10.2.2 Function Documentation	94
1	1 Class Documentation	96
	11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference	96
	11.1.1 Member Data Documentation	96
	11.2 BatteryState Class Reference	96
	11.2.1 Property Documentation	96
	11.3 BesselFilterBandPassNet Class Reference	97

11.3.1 Constructor & Destructor Documentation
11.4 BesselFilterHighPassNet Class Reference
11.4.1 Constructor & Destructor Documentation
11.5 BesselFilterLowPassNet Class Reference
11.5.1 Constructor & Destructor Documentation
11.6 ButterworthFilterBandPassNet Class Reference
11.6.1 Constructor & Destructor Documentation
11.7 ButterworthFilterHighPassNet Class Reference
11.7.1 Constructor & Destructor Documentation
11.8 ButterworthFilterLowPassNet Class Reference
11.8.1 Constructor & Destructor Documentation
11.9 CChannelTestDeviceNet Class Reference
11.9.1 Constructor & Destructor Documentation
11.9.2 Member Function Documentation
11.10 CCMOSMea_FunctionNet Class Reference
11.10.1 Constructor & Destructor Documentation
11.10.2 Member Function Documentation
11.11 CCMOSMeaDeviceNet Class Reference
11.11.1 Constructor & Destructor Documentation
11.11.2 Member Function Documentation
11.11.3 Property Documentation
11.12 CCreateFilterNet Class Reference
11.12.1 Constructor & Destructor Documentation
11.12.2 Member Function Documentation
11.12.3 Property Documentation
11.13 CDacCalibrationFunctionNet Class Reference
11.13.1 Detailed Description
11.13.2 Constructor & Destructor Documentation
11.13.3 Member Function Documentation
11.14 CDacqGroupChannelGenericSelectionNet Class Reference
11.14.1 Constructor & Destructor Documentation
11.15 CDacqGroupChannelSelectionNet Class Reference
11.15.1 Constructor & Destructor Documentation
$11.16 CDacqGroupChannelSelectionTemplateNet < DacqGroupChannelEnumTemplateNet , Dacq \\ GroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class Template Reference \ \ 12 \\ CDeviceGroupChannelInfoTemplateNet > Class TemplateNet > Class $
11.16.1 Constructor & Destructor Documentation
11.16.2 Member Function Documentation
11.17 CDeviceGroupChannelInfoGenericNet Class Reference
11.17.1 Constructor & Destructor Documentation
11.18 CDeviceGroupChannelInfoMEA2100_256Net Class Reference
11.18.1 Constructor & Destructor Documentation
11 19 CDeviceGroupChannelInfoNet Class Reference

11.19.1 Constructor & Destructor Documentation	125
11.20 CDeviceGroupChannelInfoSCUNet Class Reference	126
11.20.1 Constructor & Destructor Documentation	126
${\tt 11.21\ CDeviceGroupChannelInfoTemplateNet} < \ {\tt DacqGroupChannelEnumTemplateNet} \ > \ {\tt Class\ TemplateNet} < \ {\tt Class\ TemplateNet} > \ {\tt Class$	
plate Reference	
11.21.1 Constructor & Destructor Documentation	
11.21.2 Member Data Documentation	
11.22 CDeviceGroupChannelInfoW2100Net Class Reference	127
11.22.1 Constructor & Destructor Documentation	127
11.23 CDigOutStimulatorFunctionNet Class Reference	128
11.23.1 Detailed Description	128
11.23.2 Constructor & Destructor Documentation	129
11.23.3 Member Function Documentation	
11.24 CEncapsulatorDeviceNet Class Reference	132
11.24.1 Detailed Description	132
11.24.2 Constructor & Destructor Documentation	132
11.24.3 Member Function Documentation	133
11.25 CExternDTesterDeviceNet Class Reference	133
11.25.1 Detailed Description	133
11.25.2 Constructor & Destructor Documentation	133
11.25.3 Member Function Documentation	134
11.26 CFilterCoefficientsNet Class Reference	135
11.26.1 Constructor & Destructor Documentation	135
11.26.2 Member Function Documentation	136
11.26.3 Property Documentation	137
11.27 CFilterConfigurationNet Class Reference	137
11.27.1 Constructor & Destructor Documentation	138
11.27.2 Member Function Documentation	138
11.28 CFilterConfigurationRegisterNet Class Reference	139
11.28.1 Constructor & Destructor Documentation	140
11.28.2 Member Function Documentation	140
11.29 CFilterPropertyNet Class Reference	141
11.29.1 Constructor & Destructor Documentation	
11.29.2 Member Function Documentation	142
11.29.3 Property Documentation	142
11.30 CFluidControlDeviceNet Class Reference	
11.30.1 Detailed Description	
11.30.2 Constructor & Destructor Documentation	
11.30.3 Member Function Documentation	
11.30.4 Property Documentation	
11.31 CFYIDeviceNet Class Reference	
11.31.1 Detailed Description	

11.31.2 Constructor & Destructor Documentation
11.31.3 Property Documentation
11.32 CGenericDevelopDeviceNet Class Reference
11.32.1 Detailed Description
11.32.2 Constructor & Destructor Documentation
11.32.3 Member Function Documentation
11.33 CGilsonDeviceNet Class Reference
11.33.1 Detailed Description
11.33.2 Constructor & Destructor Documentation
11.33.3 Member Function Documentation
11.33.4 Member Data Documentation
11.34 CGrapheneASICDeviceNet Class Reference
11.34.1 Constructor & Destructor Documentation
11.34.2 Member Function Documentation
11.35 CGrapheneFunctionNet Class Reference
11.35.1 Detailed Description
11.35.2 Constructor & Destructor Documentation
11.35.3 Member Function Documentation
11.36 CHiClampDeviceNet Class Reference
11.36.1 Detailed Description
11.36.2 Constructor & Destructor Documentation
11.36.3 Property Documentation
11.37 CHLADacqNet Class Reference
11.37.1 Constructor & Destructor Documentation
11.38 CHLADeviceNet Class Reference
11.38.1 Detailed Description
11.38.2 Constructor & Destructor Documentation
11.38.3 Property Documentation
11.39 CMcsUsbDacqNet::CHWInfo Class Reference
11.39.1 Detailed Description
11.39.2 Constructor & Destructor Documentation
11.39.3 Member Function Documentation
11.40 CIntanMea_FunctionNet Class Reference
11.40.1 Constructor & Destructor Documentation
11.40.2 Member Function Documentation
11.41 CInterfaceboard2FunctionNet Class Reference
11.41.1 Detailed Description
11.41.2 Constructor & Destructor Documentation
11.41.3 Member Function Documentation
11.42 CInterfaceboardFunctionNet Class Reference
11.42.1 Detailed Description
11.42.2 Constructor & Destructor Documentation 19

11.42.3 Member Function Documentation	4
11.43 CLIH3DeviceNet Class Reference	5
11.43.1 Detailed Description	7
11.43.2 Constructor & Destructor Documentation	7
11.43.3 Member Function Documentation	7
11.43.4 Property Documentation	4
11.44 CMcsBus_AxisParametersNet Class Reference	4
11.44.1 Constructor & Destructor Documentation	5
11.44.2 Member Function Documentation	5
11.45 CMcsBus_ExtensionNet Class Reference	6
11.45.1 Constructor & Destructor Documentation	6
11.45.2 Member Function Documentation	7
11.46 CMcsBus_FYIExtensionNet Class Reference	7
11.46.1 Constructor & Destructor Documentation	7
11.46.2 Member Function Documentation	8
11.47 CMcsBus_MotorControlNet Class Reference	9
11.47.1 Constructor & Destructor Documentation	2
11.47.2 Member Function Documentation	2
11.48 CMcsBus_SensorNet Class Reference	6
11.48.1 Constructor & Destructor Documentation	8
11.48.2 Member Function Documentation	8
11.49 CMcsBus_TempSensorNet Class Reference	6
11.49.1 Constructor & Destructor Documentation	7
11.49.2 Member Function Documentation	7
11.50 CMcsBus_VoltageModeNet Class Reference	8
11.50.1 Constructor & Destructor Documentation	9
11.50.2 Member Function Documentation	0
11.51 CMcsBusNet Class Reference	3
11.51.1 Constructor & Destructor Documentation	3
11.51.2 Member Function Documentation	4
11.52 CMcsUsbDacqNet Class Reference	7
11.52.1 Detailed Description	3
11.52.2 Constructor & Destructor Documentation	3
11.52.3 Member Function Documentation	3
11.52.4 Member Data Documentation	-
11.52.5 Property Documentation	7
11.52.6 Event Documentation	7
11.53 CMcsUsbDeviceStatePushFunctionNet Class Reference	7
11.53.1 Constructor & Destructor Documentation	8
11.53.2 Member Function Documentation	8
11.53.3 Event Documentation	8
11.54 CMcsUsbDeviceStatePushNet Class Reference	8

11.54.1 Constructor & Destructor Documentation
11.54.2 Member Function Documentation
11.54.3 Event Documentation
11.55 CMcsUsbFactoryNet Class Reference
11.55.1 Constructor & Destructor Documentation
11.55.2 Member Function Documentation
11.55.3 Member Data Documentation
11.56 CMcsUsbFunctionNet Class Reference
11.56.1 Constructor & Destructor Documentation
11.56.2 Member Function Documentation
11.56.3 Member Data Documentation
11.57 CMcsUsbFunctionPointerContainer Class Reference
11.58 CMcsUsbListEntryNet Class Reference
11.58.1 Detailed Description
11.58.2 Constructor & Destructor Documentation
11.58.3 Member Function Documentation
11.58.4 Property Documentation
11.59 CMcsUsbListNet Class Reference
11.59.1 Detailed Description
11.59.2 Constructor & Destructor Documentation
11.59.3 Member Function Documentation
11.59.4 Property Documentation
11.59.5 Event Documentation
11.60 CMcsUsbNet Class Reference
11.60.1 Detailed Description
11.60.2 Constructor & Destructor Documentation
11.60.3 Member Function Documentation
11.60.4 Member Data Documentation
11.60.5 Property Documentation
11.61 CMcsUsbPointerContainer Class Reference
11.62 CMEA2100_256DacqGroupChannelSelectionNet Class Reference
11.62.1 Constructor & Destructor Documentation
11.63 CMEA2100x256FunctionNet Class Reference
11.63.1 Detailed Description
11.63.2 Constructor & Destructor Documentation
11.63.3 Member Function Documentation
11.64 CMeaAudioFunctionNet Class Reference
11.64.1 Constructor & Destructor Documentation
11.64.2 Member Function Documentation
11.65 CMeaCleanDeviceNet Class Reference
11.65.1 Detailed Description
11.65.2 Constructor & Destructor Documentation

11.65.3 Member Function Documentation
11.66 CMeaCoatDeviceNet Class Reference
11.66.1 Detailed Description
11.66.2 Constructor & Destructor Documentation
11.66.3 Member Function Documentation
11.67 CMeaDeviceNet Class Reference
11.67.1 Detailed Description
11.67.2 Constructor & Destructor Documentation
11.67.3 Member Function Documentation
11.67.4 Property Documentation
11.68 CMeaDigitalDataFunctionNet Class Reference
11.68.1 Constructor & Destructor Documentation
11.68.2 Member Function Documentation
11.69 CMeaFeedbackFunctionNet Class Reference
11.69.1 Constructor & Destructor Documentation
11.69.2 Member Function Documentation
11.70 CMealmpedanceDeviceNet Class Reference
11.70.1 Constructor & Destructor Documentation
11.70.2 Member Function Documentation
11.71 CMeasureTableDeviceNet Class Reference
11.71.1 Detailed Description
11.71.2 Constructor & Destructor Documentation
11.71.3 Property Documentation
11.72 CMeaSwitchDeviceNet Class Reference
11.72.1 Detailed Description
11.72.2 Constructor & Destructor Documentation
11.72.3 Member Function Documentation
11.73 CMeaUSBDeviceNet Class Reference
11.73.1 Detailed Description
11.73.2 Constructor & Destructor Documentation
11.73.3 Member Function Documentation
11.74 CMeFunctionNet Class Reference
11.74.1 Detailed Description
11.74.2 Constructor & Destructor Documentation
11.74.3 Member Function Documentation
11.75 CMultiBatteryChargerDeviceNet Class Reference
11.75.1 Detailed Description
11.75.2 Constructor & Destructor Documentation
11.75.3 Member Function Documentation
11.76 CMultiwellCallbackFunctionNet Class Reference
11.76.1 Detailed Description
11.76.2 Constructor & Destructor Documentation

11.76.3 Member Function Documentation	383
11.76.4 Event Documentation	383
11.77 CMultiwellDeviceNet Class Reference	384
11.77.1 Detailed Description	385
11.77.2 Constructor & Destructor Documentation	385
11.77.3 Member Function Documentation	386
11.78 CMultiwellOptoStimFunctionNet Class Reference	391
11.78.1 Detailed Description	392
11.78.2 Constructor & Destructor Documentation	392
11.78.3 Member Function Documentation	392
11.79 CNF_GenDeviceNet Class Reference	396
11.79.1 Constructor & Destructor Documentation	396
11.79.2 Member Function Documentation	397
11.80 COctoPotDeviceNet Class Reference	397
11.80.1 Constructor & Destructor Documentation	398
11.80.2 Member Function Documentation	398
11.81 COkuvisionStimulatorDeviceNet Class Reference	401
11.81.1 Constructor & Destructor Documentation	402
11.81.2 Member Function Documentation	402
11.82 CPatchServerDeviceNet Class Reference	404
11.82.1 Detailed Description	405
11.82.2 Constructor & Destructor Documentation	405
11.82.3 Property Documentation	405
11.83 CPathIdentDeviceNet Class Reference	405
11.83.1 Constructor & Destructor Documentation	406
11.83.2 Member Function Documentation	406
11.84 CPedoterDeviceNet Class Reference	406
11.84.1 Detailed Description	407
11.84.2 Constructor & Destructor Documentation	407
11.84.3 Member Function Documentation	407
11.85 CPeristalticPumpDeviceNet Class Reference	408
11.85.1 Detailed Description	408
11.85.2 Constructor & Destructor Documentation	408
11.85.3 Property Documentation	409
11.86 CPgaDeviceNet Class Reference	409
11.86.1 Constructor & Destructor Documentation	409
11.86.2 Member Function Documentation	410
11.87 CPositionIIDeviceNet Class Reference	411
11.87.1 Detailed Description	413
11.87.2 Constructor & Destructor Documentation	413
11.87.3 Member Function Documentation	413
11.87.4 Property Documentation	420

11.88 CPositionImpDeviceNet Class Reference	0
11.88.1 Detailed Description	0
11.88.2 Constructor & Destructor Documentation	1
11.88.3 Member Function Documentation	1
11.89 CPPCDeviceNet Class Reference	3
11.89.1 Constructor & Destructor Documentation	3
11.89.2 Property Documentation	4
11.90 CPPCFunctionNet Class Reference	4
11.90.1 Detailed Description	5
11.90.2 Constructor & Destructor Documentation	5
11.90.3 Member Function Documentation	6
11.91 CPPS_DeviceNet Class Reference	3
11.91.1 Constructor & Destructor Documentation	3
11.91.2 Property Documentation	3
11.92 CPPS_FunctionNet Class Reference	4
11.92.1 Constructor & Destructor Documentation	4
11.92.2 Member Function Documentation	5
11.93 CPPSDeviceNet Class Reference	8
11.93.1 Detailed Description	8
11.93.2 Constructor & Destructor Documentation	9
11.94 CProgramPressureCurveNet Class Reference	9
11.94.1 Detailed Description	9
11.94.2 Constructor & Destructor Documentation	9
11.94.3 Member Function Documentation	0
11.95 CPulseGeneratorFunctionNet Class Reference	0
11.95.1 Detailed Description	1
11.95.2 Constructor & Destructor Documentation	1
11.95.3 Member Function Documentation	2
11.96 CRadioControledDevicesNet Class Reference	3
11.96.1 Constructor & Destructor Documentation	4
11.96.2 Member Function Documentation	4
11.97 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference	5
11.97.1 Constructor & Destructor Documentation	5
11.97.2 Member Function Documentation	5
11.97.3 Member Data Documentation	6
11.98 CRetinaLedDeviceNet Class Reference	6
11.98.1 Constructor & Destructor Documentation	7
11.98.2 Member Function Documentation	7
11.99 CRFFunctionNet Class Reference	8
11.99.1 Detailed Description	9
11.99.2 Constructor & Destructor Documentation	9
11.99.3 Member Function Documentation 44	9

11.100 CRobo_FYIProgram_FunctionNet Class Reference
11.100.1 Constructor & Destructor Documentation
11.100.2 Member Function Documentation
11.101 CRobo_FYITemp_FunctionNet Class Reference
11.101.1 Constructor & Destructor Documentation
11.101.2 Member Function Documentation
11.102 CRoboDacqNet Class Reference
11.102.1 Constructor & Destructor Documentation
11.102.2 Member Function Documentation
11.102.3 Member Data Documentation
11.103 CRoboDeviceNet Class Reference
11.103.1 Detailed Description
11.103.2 Constructor & Destructor Documentation
11.103.3 Member Function Documentation
11.103.4 Member Data Documentation
11.103.5 Property Documentation
11.103.6 Event Documentation
11.104 CRoboFluidDeviceNet Class Reference
11.104.1 Constructor & Destructor Documentation
11.104.2 Member Function Documentation
11.104.3 Member Data Documentation
11.104.4 Property Documentation
11.105 CRoboInjectDeviceNet Class Reference
11.105.1 Detailed Description
11.105.2 Constructor & Destructor Documentation
11.106 CRoboocyte2DeviceNet Class Reference
11.106.1 Detailed Description
11.106.2 Constructor & Destructor Documentation
11.106.3 Member Function Documentation
11.107 CRoboStatorDeviceNet Class Reference
11.107.1 Constructor & Destructor Documentation
11.107.2 Member Function Documentation
11.107.3 Property Documentation
11.108 CSafeISDeviceNet Class Reference
11.108.1 Detailed Description
11.108.2 Constructor & Destructor Documentation
11.108.3 Member Function Documentation
11.108.4 Property Documentation
11.109 CSCUDacqGroupChannelSelectionNet Class Reference
11.109.1 Constructor & Destructor Documentation
11.110 CSCUFunctionNet Class Reference
11.110.1 Detailed Description 504

11.110.2 Constructor & Destructor Documentation	5
11.110.3 Member Function Documentation	5
11.110.4 Event Documentation	9
11.111 CSerialPortNet Class Reference	9
11.111.1 Constructor & Destructor Documentation	9
11.111.2 Member Function Documentation	0
11.112 CStg200xBasicNet Class Reference	0
11.112.1 Detailed Description	:5
11.112.2 Constructor & Destructor Documentation	:5
11.112.3 Member Function Documentation	6
11.113 CStg200xDownloadBasicNet Class Reference	9
11.113.1 Detailed Description	0
11.113.2 Member Function Documentation	1
11.113.3 Property Documentation	6
11.114 CStg200xDownloadNet Class Reference	6
11.114.1 Detailed Description	7
11.114.2 Constructor & Destructor Documentation	7
11.114.3 Member Function Documentation	8
11.114.4 Event Documentation	3
11.115 CStimulusFunctionNet Class Reference	3
11.115.1 Constructor & Destructor Documentation	5
11.115.2 Member Function Documentation	5
11.115.3 Event Documentation	4
11.116 CSw2to64DeviceNet Class Reference	4
11.116.1 Detailed Description	4
11.116.2 Constructor & Destructor Documentation	5
11.116.3 Member Function Documentation	5
11.117 CTcxDeviceNet Class Reference	6
11.117.1 Detailed Description	8
11.117.2 Constructor & Destructor Documentation	8
11.117.3 Member Function Documentation	9
11.118 CTEERFunctionNet Class Reference	9
11.118.1 Detailed Description	1
11.118.2 Constructor & Destructor Documentation	1
11.118.3 Member Function Documentation	2
11.119 CTEERMachineDeviceNet Class Reference	0
11.119.1 Constructor & Destructor Documentation	
11.119.2 Property Documentation	1
11.120 CUsbDeviceConfigurationFunctionNet Class Reference	
11.120.1 Detailed Description	1
11.120.2 Constructor & Destructor Documentation	1
11.120.3 Member Function Documentation 61	2

11.121 CUsbExceptionNet Class Reference
11.121.1 Detailed Description
11.121.2 Constructor & Destructor Documentation
11.121.3 Property Documentation
11.122 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference 613
11.122.1 Constructor & Destructor Documentation
11.122.2 Member Data Documentation
11.123 CW2100_FunctionNet Class Reference
11.123.1 Constructor & Destructor Documentation
11.123.2 Member Function Documentation
11.123.3 Property Documentation
11.124 CW2100_StimulatorFunctionNet Class Reference
11.124.1 Constructor & Destructor Documentation
11.124.2 Member Function Documentation
11.124.3 Member Data Documentation
11.124.4 Event Documentation
11.125 CW2100DacqGroupChannelSelectionNet Class Reference
11.125.1 Constructor & Destructor Documentation
11.126 CWarnerUssingDeviceNet Class Reference
11.126.1 Detailed Description
11.126.2 Constructor & Destructor Documentation
11.126.3 Property Documentation
11.127 CWarnerUssingFunctionNet Class Reference
11.127.1 Detailed Description
11.127.2 Constructor & Destructor Documentation
11.127.3 Member Function Documentation
11.128 CWarnerValveControllerDeviceNet Class Reference
11.128.1 Detailed Description
11.128.2 Constructor & Destructor Documentation
11.128.3 Member Function Documentation
11.128.4 Event Documentation
11.129 CWarnerValveControllerDeviceTesterFunctionNet Class Reference 67
11.129.1 Detailed Description
11.129.2 Constructor & Destructor Documentation
11.129.3 Member Function Documentation
11.130 CWClassicFunctionNet Class Reference
11.130.1 Constructor & Destructor Documentation
11.130.2 Member Function Documentation
11.131 CWirelessBaseFunctionNet Class Reference
11.131.1 Constructor & Destructor Documentation
11.131.2 Member Function Documentation
11 132 Device IdNat Struct Reference

11.132.1 Detailed Description
11.132.2 Constructor & Destructor Documentation
11.132.3 Member Function Documentation
11.132.4 Member Data Documentation
11.133 DigitalSource< digitalsourceenum > Class Template Reference
11.133.1 Constructor & Destructor Documentation
11.133.2 Member Function Documentation
11.133.3 Property Documentation
11.134 DigitalSourceGeneral Class Reference
11.134.1 Constructor & Destructor Documentation
11.134.2 Member Function Documentation
11.134.3 Property Documentation
11.135 DriverVersionNet Class Reference
11.135.1 Detailed Description
11.135.2 Constructor & Destructor Documentation
11.135.3 Member Function Documentation
11.136 FirmwareDestinationNames Class Reference
11.136.1 Member Data Documentation
11.137 HeadStageIDType Class Reference
11.137.1 Member Enumeration Documentation
11.137.2 Constructor & Destructor Documentation
11.137.3 Member Function Documentation
11.137.4 Property Documentation
11.138 HeadstageIDTypeObject Class Reference
11.138.1 Constructor & Destructor Documentation
11.138.2 Member Function Documentation
11.138.3 Member Data Documentation
11.138.4 Property Documentation
11.139 HeadStageIDTypeState Class Reference
11.139.1 Property Documentation
11.140 mkfilterNet Class Reference
11.140.1 Member Function Documentation
11.141 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference
11.141.1 Member Function Documentation
11.142 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference 709
11.142.1 Member Function Documentation
11.143 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference
11.143.1 Constructor & Destructor Documentation
11.143.2 Member Function Documentation
11.143.3 Member Data Documentation
11.144 CMeaAudioFunctionNet::s_setaudionet Struct Reference
11 1// 1 Member Data Documentation 71

11.145 CStimulusFunctionNet::SidebandData Class Reference	711
11.145.1 Constructor & Destructor Documentation	711
11.145.2 Property Documentation	711
11.146 StgStatusNet Class Reference	712
11.146.1 Member Function Documentation	712
11.146.2 Member Data Documentation	712
11.147 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference	712
11.147.1 Constructor & Destructor Documentation	713
11.147.2 Property Documentation	713
11.148 CMcsUsbFactoryNet::TargetSelectionInfoNet Class Reference	714
11.148.1 Member Data Documentation	714
11.149 usbSetupPacket_t Class Reference	714
11.149.1 Member Data Documentation	714
11.150 W2100_StimulusParametersNet Struct Reference	715
11.150.1 Member Data Documentation	715
Index	717

1 McsUsbNet.dll for MCS USB devices

1.1 Introduction

This DLL provides the .NET interface to MCS devices

The most important options are accessing our stimulator and data acquisition devices:

- STG200x & STG400x STimulus Generator
- · Data ACQuisition (DACQ) Devices

See here for a list of our other devices: Device Classes.

And here for a list of function classes addressing groups of features that might be shared between different devices: Function Classes.

1.2 System requirements

The DLL can be used with any .NET compatible language.

The DLL needs the .NET Framework 4.7.2.

It requires the Microsoft Visual C++ Redistributable for Visual Studio 2019 to be installed.

It also requires the **USB driver** to be installed.

The simplest way to achieve this is to install the latest **Multi Channel Experimenter** setup (will install 64bit redistributable).

All examples assume that the Mcs.Usb namespace is loaded:

sing namespace Mcs.Usb;

Include the file McsUsbNet.dll into the references of your project.

1.3 Connecting to an MCS device

A connection to a DAQ device is established by Mcs.Usb.CMcsUsbNet.Connect. When this function is called without argument, the first DAQ device found on the USB bus is used:

```
CMcsUsbNet device = new CMcsUsbNet();
device.Connect();
```

When more than one DAQ device of the specific type is connected, you can use the Mcs.Usb.CMcsUsbListNet class to get a list of available devices:

```
CMcsUsbListNet usblist = new CMcsUsbListNet(DeviceEnumNet.MCS_DEVICE_USB);
var entry = usblist.GetUsbListEntry((uint)0);
CMcsUsbNet device = new CMcsUsbNet();
device.Connect(entry);
```

After you are finished with the device, you can disconnect the device object from the device by: device.Disonnect();

2 Device Classes

- For FluidControl device see MCS FluidControl
- For SW2TO64 device see MCS-USB-Sw2to64
- For TCx device see Mcs.Usb.CTcxDeviceNet

2.1 The MCS FluidControl Device

2.1.1 Introduction

The FluidControl Device can control up to 24 valves. The nominal voltage is 24V.

8 TTL level digital output ports are available and 8 TTL inputs can be read in.

The device has 8 ADC inputs with a rage from 0V to 3.3V.

2.1.2 Access to the FluidControl device

status = m_fluidcontrol->Connect();

For connecting to a FluidControl device see Connecting to an MCS device.*

CFluidControlDevice* m_dacq;

m_fluidcontrol = new CFluidControlDevice;

The valves are controlled with the CFluidControlDevice::SetValve call. The argument given is a bit pattern of all valves which should be open.

The digital outputs can be controlled with the CFluidControlDevice::SetDigout call. Again, a bit pattern of all digital output pins which should be set to a logic high level is given as an argument.

The current state of the valves and the digital outputs can be read back with the CFluidControlDevice::GetValve and CFluidControlDevice::GetDigout

The command to read an ADC-Channel is CFluidControlDevice::GetAdc. Here the channelnummer which should be read in is given as an argument and the return value is the current Adc level.

The state of the digital inputs is read with the CFluidControlDevice::GetDigin call. Here the return value is the bit pattern of the digital inputs.

The connection to the device is closed with the CFluidControlDevice::Disconnect call.

2.2 MCS-USB-Sw2to64 device

The class Mcs.Usb.CSw2to64DeviceNet controls the setting of the switches in the MCS-USB-Sw2to64 device.

First construct an object of the class:

```
CSw2to64DeviceNet device = new CSw2to64DeviceNet();
```

For connecting to an MCS-USB-Sw2to64 device see Connecting to an MCS device.

To get the number of channels the device handles:

```
int number = device.GetNumber()
```

Set all channel switches at once:

```
byte z = 1;
byte[] pattern = new byte[number];
for(int i = 0;i < number;i++)
{
   pattern[i] = z; // pattern you want to switch this channel to
}
device.SetChannels(pattern);</pre>
```

Get all channel switches at once:

```
byte[] pattern = device.GetChannels();
```

Set one channel switch:

```
ushort index = 10;
byte pattern = 1;
device.SetChannel(index, pattern)
```

Get one channel switch:

```
ushort index = 10;
byte pattern = device.GetChannel(index);
```

3 Function Classes

- Mcs.Usb.CCMOSMea_FunctionNet
- Mcs.Usb.CDacCalibrationFunctionNet
- Mcs.Usb.CDigOutStimulatorFunctionNet
- Mcs.Usb.CGrapheneFunctionNet
- Mcs.Usb.CIntanMea_FunctionNet
- Mcs.Usb.CInterfaceboard2FunctionNet
- Mcs.Usb.CInterfaceboardFunctionNet
- Mcs.Usb.CMcsBus_MotorControlNet
- Mcs.Usb.CMcsBus_VoltageModeNet
- Mcs.Usb.CMcsBus_AxisParametersNet
- Mcs.Usb.CMcsBus_SensorNet
- Mcs.Usb.CMcsBus_TempSensorNet
- Mcs.Usb.CMcsBus_ExtensionNet
- Mcs.Usb.CMcsBus_FYIExtensionNet
- Mcs.Usb.CMcsUsbDeviceStatePushFunctionNet
- Mcs.Usb.CMEA2100x256FunctionNet

- Mcs.Usb.CMeaAudioFunctionNet
- Mcs.Usb.CMeaDigitalDataFunctionNet
- Mcs.Usb.CMeaFeedbackFunctionNet
- · Mcs.Usb.CMeFunctionNet
- Mcs.Usb.CMultiwellCallbackFunctionNet
- Mcs.Usb.CMultiwellOptoStimFunctionNet
- Mcs.Usb.CPPCFunctionNet
- Mcs.Usb.CPPS_FunctionNet
- · Mcs.Usb.CPulseGeneratorFunctionNet
- Mcs.Usb.CRFFunctionNet
- Mcs.Usb.CRobo_FYITemp_FunctionNet
- Mcs.Usb.CRobo_FYIProgram_FunctionNet
- Mcs.Usb.CSCUFunctionNet
- Mcs.Usb.CStimulusFunctionNet
- Mcs.Usb.CTEERFunctionNet
- Mcs.Usb.CW2100 FunctionNet
- Mcs.Usb.CW2100_StimulatorFunctionNet
- Mcs.Usb.CWarnerUssingFunctionNet
- Mcs.Usb.CWarnerValveControllerDeviceTesterFunctionNet
- Mcs.Usb.CWClassicFunctionNet
- Mcs.Usb.CWirelessBaseFunctionNet

4 Data ACQuisition (DACQ) Devices

There are different device types of (MEA) data acquisition (DACQ) devices. All of them are supported by this class.

This library does **not** support the writing of the MCD (MC_Rack), MSRD (Multi Channel Experimenter) or HDF5 file format!

The class Mcs.Usb.CMeaDeviceNet is the base class for DACQ devices.

The base class Mcs.Usb.CMeaDeviceNet constructs actually the underlying classes for USB-MEA devices (Mcs.Usb.CMeaUSBDeviceNet).

```
.
CMeaDeviceNet device = new CMeaDeviceNet(McsBusTypeEnumNet.MCS_USB_BUS, OnChannelData, OnError);
```

For connecting to a DACQ device see Connecting to an MCS device.

Get the number of available analog hardware channels and set the number of channels to the maximum.

```
int hwchannels;
device.HWInfo().GetNumberOfHWADCChannels(out hwchannels);
device.SetNumberOfChannels(hwchannels);
int samplingrate = 1000;
device.SetSamplerate(samplingrate, 1, 0);
device.EnableDigitalIn(true, 0);
```

5 The MCS Robo Device 5

Get the layout to know how the data look like that you receive

```
int ana, digi, che, tim, block;
device.GetChannelLayout(out ana, out digi, out che, out tim, out block);
```

For the Mcs.Usb.OnChannelData callback function you have to provide a definition of the channels you want to receive.

```
bool[] selChannels = new bool[block];
for (int i = 0; i < block; i++)
{
    selChannels[i] = true; // With true channel i is selected
    // selChannels[i] = false; // With false the channel i is deselected
}
channelblocksize = samplingrate / 10;
// queue size and threshold should be selected carefully
device.SetSelectedChannels(selChannels, 10 * channelblocksize, channelblocksize);</pre>
```

The Mcs.Usb.OnChannelData callback function gets a callback for each channelblock that is defined. In this example a callback for each channel.

```
void OnChannelData(CMcsUsbDacqNet d, int cbHandle, int numSamples)
{
   int size_ret;
   ushort[] channeldata = device.ChannelBlock_ReadFramesUI16(CbHandle, numSamples, out size_ret);
}
void OnError(String msg, int info)
{
   MessageBox.Show("Mea Device Error: " + msg);
}
```

see MEA Recording in the Examples directory.

5 The MCS Robo Device

5.1 Introduction

Up to now two MCS devices exist that base on the Robo platform.

- The MCS Roboinject device is controlled by the Mcs.Usb.CRobolnjectDeviceNet class.
- The MCS Roboocyte2 device is controlled by the Mcs.Usb.CRoboocyte2DeviceNet class.

Both classes are derived from Mcs.Usb.CRoboDeviceNet

6 STG200x & STG400x STimulus Generator

6.1 Introduction

The STG200x & STG400x Series Stimulus Generators have two distinct modes of operation, the Download mode and the Streaming mode.

6.2 Download mode

The Download mode is the "classic" mode of operation, as used by the MC Stimulus software. In this mode, one or multiple waveforms are defined in PC memory and downloaded to the STG. The waveforms are stored in STG device onboard memory and can be sent to the analog and sync outputs once or multiple times. The STG can operate independently from the PC (without computer connection) after the download. Output is triggered either by the front panel start/stop button, the digital trigger inputs or under software control.

In the Download mode, there are up to eight independent triggers available (depending on the device). The user can assign each of the analog outputs and sync (digital) outputs to any of the triggers.

The analog output waveform is stored sample by sample in the STG memory. To reduce memory usage, this data can be compressed: whenever a given output value is to be held for more than one sample period, it has only to be given once. The user can define the number of sample periods for that a pattern should remain active. Compression is done for each channel independently of the others, thus the algorithm to compress the data is very easy to implement.

A new feature of the Download mode is the segmentation of the STG memory. The onboard memory can be devided into up to 100 segments. Each segment can hold its own waveform pattern. Under software control, the user can switch between the defined segments within milliseconds. Another option is to use the four trigger inputs to select between four predefined segments. This option is accessible from the MC_Stimulus Software as the "Multi-File mode", and can start each of up to four defined waveforms within microseconds. This feature allows a predefinied flexible response (feedback) to recorded data.

Mcs.Usb.CStg200xDownloadNet is the class for using the STG in download mode.

6.2.1 Memory Layout and Trigger Setup

The class to be used for the Download mode is Mcs.Usb.CStg200xDownloadNet, which is derived from Mcs.Usb.CStg200xBasicNet. You can add a poll handler delegate (Mcs.Usb.OnStg200xPollStatus) to the constructor Mcs.Usb.CStg200xDownloadNet.

For connecting to an STG see Connecting to an MCS device.

To use the Download mode, the memory layout of the STG200x can be set up, if the default is not sufficient. The total amount of memory available in the STG is obtained by the Mcs.Usb.CStg200xDownloadNet.GetTotalMemory call. With Mcs.Usb.CStg200xDownloadNet.SegmentDefine the segment sizes are assigned.

Next, for each segment, one has to assign the amount of memory to be used for each channel and sync output. This is done by Mcs.Usb.CStg200xDownloadBasicNet.SetCapacity. Its arguments contain a list of memory sizes, with one entry per channel and one entry per sync output. Again, the total memory assigned to the channels and sync outputs must not exceed the memory assigned to the segment.

6.2 Download mode 7

}

Before the STG can start, the trigger has to be configured. This is done by the Mcs.Usb.CStg200xDownloadNet.SetupTrigger call. Its arguments are a list of channelmaps, syncoutmaps and repeats, one for each of the four available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and syncout 1 to trigger 1 and channel 3 to trigger 2 use:

```
uint32_t TriggerInputs = device.GetNumberOfTriggerInputs();
uint[] channelmap = new uint[TriggerInputs];
uint[] syncoutmap = new uint[TriggerInputs];
uint[] repeat = new uint[TriggerInputs];
for (int i = 0; i < TriggerInputs; i++)
{
    channelmap[i] = 0;
    syncoutmap[i] = 0;
    repeat[i] = 0;
}
// Trigger 0
    channelmap[0] = 1; // Channel 1
    syncoutmap[0] = 1; // Syncout 1
    repeat[0] = 0; // forever
// Trigger 1
    channelmap[1] = 4; // Channel 3
    device.SetupTrigger(channelmap, syncoutmap, repeat);</pre>
```

For the STG400x series you have to set the output mode of the channels. Mcs.Usb.CStg200xDownloadNet.SetVoltageMode interprets the values as voltages. Mcs.Usb.CStg200xDownloadNet.SetCurrentMode as currents.

```
// Only meaningfull for STG400x
device.SetVoltageMode();
```

For each segment, data can be sent to each of the defined channels and sync outputs using the Mcs.Usb.CStg200xDownloadNet.SendChannelData and Mcs.Usb.CStg200xDownloadNet.SendSyncData calls. channeldata and syncdata are a list of analog and digital samples as a list of two byte values (unsigned short). Multiple calls to Mcs.Usb.CStg200xDownloadNet.SendSyncData to the same channel append data to that channel.

If the Multi-File mode of the STG is enabled using the Mcs.Usb.CStg200xDownloadNet.EnableMultiFileMode call, the four trigger inputs are used to switch between four segments. A hardware trigger signal (TTL) on trigger input 1 selects the first segment and starts all pulses in this segment. Thus with the Multi-File mode, one can predefine four stimulus patterns and switch between them without a connection to the PC.

The STG200x series has an analog resolution of 13 bits, thus the analog data contains the information in bits 0 to 12 of each sample. Bits 13 to 15 have to be 0.

```
int DACResolution = device.GetDACResolution();
// Data for Channel 0
    device.ClearChannelData(0);
    double factor = 0.1;
    const int 1 = 1000;
   ushort[] pData = new ushort[1];
        Uint64_t[] tData = new Uint64_t[1];
        for (int i = 0; i < 1; i++)
             // calculate Sin-Wave
            double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
                Math.Sin(2.0 * (double)i * Math.PI / (double)1);
             // calculate sign
            pData[i] = sin >= 0 ? (ushort)sin : (ushort)((int)Math.Abs(sin) +
                 (int)Math.Pow(2, DACResolution - 1));
            tData[i] = (Uint64_t)20; // duration in \mus
        device. SendChannelData(0, pData, tData);
// Data for Channel 3
    device.ClearChannelData(2);
    double factor = 0.1;
    const int 1 = 700;
    // without compression
    ushort[] pData = new ushort[1];
    Uint64_t[] tData = new Uint64_t[1];
    for (int i = 0; i < 1; i++)
        // calculate Sin-Wave
        double sin = factor * (Math.Pow(2, DACResolution - 1) - 1.0) *
           Math.Sin(2.0 * (double)i * Math.PI / (double)1);
```

Start the trigger by pushing the front button or by software

```
// Start Trigger 1 and 2
device.SendStart(1 + 2); // Trigger 1 und 2
```

see the StgDownloadExampleNet in the example directory.

6.3 Streaming mode

The other mode of operation is the Streaming mode. Here the analog output is sent to the STG device in "real time". The PC has to be connected to the STG all the time. The data that is sent to the analog output is downloaded from the PC to the STG on the fly.

The Streaming mode is useful for applications where flexible feedback is needed as well for applications where very long waveforms which are not repeated (such as white noise) are used.

The Streaming mode works by use of two ring buffers which hold data. One is in PC memory and managed by the DLL, and one is in on-board STG memory. Data is transferred from PC memory to the STG via the USB bus in time slices of one millisecond.

The user can define both the size of the ring buffer in DLL memory and in the STG memory. Once the Streaming mode is started, the STG request data from the PC. The data rate from PC to STG is variable and controlled by the STG. The STG request data from the PC at a rate to keep its internal ringbuffer at about half full.

It is the responsibility of the user to keep the ring buffer in the memory of the PC filled, so the DLL can supply sufficient data to the STG. To do so, the Windows DLL allows to define a "callback" function which is called whenever new data is needed, or more precise, as soon as the ring buffer in the memory of the PC falls below the user defined threshold.

Small buffers have the advantage of a low latency between data generation in the callback funtion and its output as a analog signal from the STG. However for low latency to work, the user-written callback function has to be fast and to produce a steady flow of data.

In the Streaming mode, all triggers are available as well. Each of the eight analog and sync outputs can be assigned to one of the triggers.

The output rate is user defined with a maximum of 50 kHz

Mcs.Usb.CStg200xStreamingNet is the class for using the STG in streaming mode.

6.3 Streaming mode 9

6.3.1 Memory Layout and Trigger Setup

With the constructor for Mcs.Usb.CStg200xStreamingNet.CStg200xStreamingNet, the name of the callback function for the data handler is provided. The data handler function is called automatically, whenever the STG needs new data. This data is first written to a ring buffer in the memory of the PC. The size for this ring buffer is defined as first argument in the constructor. The user provided delegate gets the trigger number which needs new data as argument

```
CStg200xStreamingNet device = new CStg200xStreamingNet(10000, dataHandler, errorHandler);
```

The callback funtion, which is defined in the constructor, is called whenever the STG needs new data for a trigger, or more precise, whenever the ring buffer in PC memory falls below the defined threshold.

The user can query the amount of space available for queuing by use of the Mcs.Usb.CStg200xStreamingNet. ← GetDataQueueSpace call. Its return value is the number of samples that can be send to the STG.

User code is required to fill an array analog and sync out data, sample by sample for up to the maximum number of samples as obtained by Mcs.Usb.CStg200xStreamingNet.GetDataQueueSpace or Mcs.Usb.CStg200xStreaming

Net.GetSyncoutQueueSpace.

The values for the analog outputs are 16 bits signed integers. The lower bits are trunctated according to the resolution of the STG. This behaviour is different to the behaviour in download mode.

Note: Compression as described in the download mode can NOT be used for the streaming mode.

The new data is sent to the STG by using the Mcs.Usb.CStg200xStreamingNet.EnqueueData call.

```
void dataHandler(uint32_t trigger)
    double factor = 1;
if (trigger == 0) // Callback for Trigger 1
          {// Handle Channel 1
               uint32_t channel = 0;
                    uint32 t space = device.GetDataOueueSpace(channel);
                    if (space < 1000)
                    short[] data = new short[1000];
                    for (int i = 0; i < 1000; i++)
                         // Calc Sin-Wave (16 bits) lower bits will be removed according resolution double sin = factor \star (Math.Pow(2, 16 - 1) - 1.0) \star Math.Sin(2.0 \star (double)i \star Math.PI / (double)1000);
                         data[i] = (short)sin;
                    uint32_t engueued = device.EngueueData(channel, data);
          {// Handle Channel 3
               uint32_t channel = 2;
               for (; ; )
                    uint32_t space = device.GetDataQueueSpace(channel);
                    if (space < 700)</pre>
                         break:
                    short[] data = new short[700];
                    for (int i = 0; i < 700; i++)
                         // Calc Sin-Wave (16 bits) lower bits will be removed according resolution double sin = factor * (Math.Pow(2, 16 - 1) - 1.0) * Math.Sin(2.0 * (double)i * Math.PI / (double)700);
                         data[i] = (short)sin;
                    uint32_t enqueued = device.EnqueueData(channel, data);
          {// Handle Syncout 1
               uint32_t channel = 0;
               for (; ; )
                    uint32_t space = device.GetSyncoutQueueSpace(channel);
                    if (space < 1000)</pre>
                         break:
                    ushort[] data = new ushort[1000];
                    for (int i = 0; i < 1000; i++)</pre>
```

For connecting to an STG device see Connecting to an MCS device.

With enabling or disabling the continuous mode it can be selected how the STG handles an "out of data" situation.

When Mcs.Usb.CStg200xStreamingNet.EnableContinousMode is used, the STG does not stop when it runs out of data, but it keeps running and sends a zero voltage to its outputs.

When Mcs.Usb.CStg200xStreamingNet.DisableContinousMode is used, the STG stops when it runs out of data. It has to be retriggered to resume the output.

```
device.EnableContinousMode();
```

Mcs.Usb.CStg200xStreamingNet.SetOutputRate is used to set the sampling rate. device.SetOutputRate(50000);

To use the Streaming mode, the memory layout of the STG has to be set up. To total amount of memory available in the STG is obtained by the Mcs.Usb.CStg200xStreamingNet.GetTotalMemory call.

This memory can be assigned to four ring buffers (one per trigger) which buffer the data received from the PC via USB cable. This is done with the CStg200xStreaming::SetCapacity call. The total amount of memory must not exceed the total memory size as obtained by Mcs.Usb.CStg200xStreamingNet.GetTotalMemory.

This internal ring buffer is crucial for proper operation of the Streaming mode. The size of the ring buffer determines the latency of the Streaming mode. The firmware of the STG requests data from the PC in order to keep the ring buffer about half full. Thus the average latency is:

```
latency = (ringbuffersize in bytes/4) / output rate
```

If the ring buffer size is too big, the latency of the STG might be too long. If the ring buffer size is too low, an overflow or underflow of data in the STG ringbuffer might occur, resulting in data jumps of the output signals or the "out of data" situation described erlier.

The following example divides the total memory equally amoung the four triggers:

Before the STG can start, the trigger has to be configured. This is done by the Mcs.Usb.CStg200xStreaming
Net.SetupTrigger call. Its arguments are a list of channelmaps, syncoutmaps, digoutmap, autostart and callback
_threshold, with one entry for each of the available triggers. channelmap is a bitmap, each bit representing one of the available channels. To assign channel 1 and 3 and syncout 1 to trigger 1 use:

```
uint32_t ntrigger = device.GetNumberOfTriggerInputs();  // obtain number of triggers in this STG
uint[] channelmap = new uint[ntrigger];
uint[] syncoutmap = new uint[ntrigger];
uint[] digoutmap = new uint[ntrigger];
uint[] autostart = new uint[ntrigger];
uint[] callback_threshold = new uint[ntrigger];
for (int i = 0; i < ntrigger; i++)
{
    channelmap[i] = 0;</pre>
```

7 Namespace Index 11

```
syncoutmap[i] = 0;
digoutmap[i] = 0;
autostart[i] = 0;
autostart[i] = 0;
callback_threshold[i] = 0;
}
channelmap[0] = 0x1 + 0x4; // Channel 1 und Channel 3 to Trigger 1
syncoutmap[0] = 0x1; // Syncout 1 to Trigger 1
autostart[0] = 1;
callback_threshold[0] = 50; // 50% of buffer size
device.SetupTrigger(channelmap, syncoutmap, digoutmap, autostart, callback_threshold);
device.StartLoop();
System.Threading.Thread.Sleep(1000); // Give StartLoop some time
```

Start Trigger by pushing the front button or by Software

device.SendStart(1);

see the StgStreamingExampleNet in the example directory.

7 Namespace Index

7.1 Namespace List

Here is a list of all namespaces with brief descriptions:

Mcs	22
Mcs::Usb	23

8 Hierarchical Index

8.1 Class Hierarchy

This inheritance list is sorted roughly, but not completely, alphabetically:

CW2100_FunctionNet::AudioChannelsNet	96
BatteryState	96
CCreateFilterNet	116
BesselFilterBandPassNet	97
BesselFilterHighPassNet	97
BesselFilterLowPassNet	98
ButterworthFilterBandPassNet	99
ButterworthFilterHighPassNet	99
ButterworthFilterLowPassNet	100
${\tt CDeviceGroupChannelInfoTemplateNet} < {\tt DacqGroupChannelEnumTemplateNet} >$	126
${\bf CDeviceGroupChannelInfoTemplateNet} < {\bf DacqGroupChannelEnumNet} >$	126
CDeviceGroupChannelInfoNet	125

${\sf CDeviceGroupChannelInfoTemplateNet} < {\sf int} >$	126
CDeviceGroupChannelInfoGenericNet	124
CDeviceGroupChannelInfoTemplateNet< MEA2100_256DacqGroupChannelEnumNet >	126
CDeviceGroupChannelInfoMEA2100_256Net	124
CDeviceGroupChannelInfoTemplateNet< SCUDacqGroupChannelEnumNet >	126
CDeviceGroupChannelInfoSCUNet	126
CDeviceGroupChannelInfoTemplateNet< W2100DacqGroupChannelEnumNet >	126
CDeviceGroupChannelInfoW2100Net	127
CFilterCoefficientsNet	135
CFilterPropertyNet	141
CMcsUsbDacqNet::CHWInfo	186
CMcsUsbFunctionNet	308
CDacqGroupChannelSelectionTemplateNet< W2100DacqGroupChannelEnumNet, W2100← DacqGroupChannelEnum, CDeviceGroupChannelInfoW2100Net >	121
CW2100DacqGroupChannelSelectionNet	628
${\bf CDacqGroupChannelSelectionTemplateNet} < {\bf DacqGroupChannelEnumNet}, {\bf DacqGroup} \leftarrow \\ {\bf ChannelEnum, CDeviceGroupChannelInfoNet} >$	121
CDacqGroupChannelSelectionNet	121
${\tt CDacqGroupChannelSelectionTemplateNet} {< \tt int, int, CDeviceGroupChannelInfoGenericNet} {< \tt int, int, int, int, int, int, int, int,$	121
CDacqGroupChannelGenericSelectionNet	120
${\bf CDacqGroupChannelSelectionTemplateNet} < {\bf SCUDacqGroupChannelEnumNet}, {\bf SCUDacq} \leftarrow {\bf GroupChannelEnum, CDeviceGroupChannelInfoSCUNet} >$	121
CSCUDacqGroupChannelSelectionNet	501
CDacqGroupChannelSelectionTemplateNet< MEA2100_256DacqGroupChannelEnumNet, MEA2100_256DacqGroupChannelEnum, CDeviceGroupChannelInfoMEA2100_256Net >	121
CMEA2100_256DacqGroupChannelSelectionNet	339
CCMOSMea_FunctionNet	102
CDacCalibrationFunctionNet	118
$\textbf{CDacqGroupChannelSelectionTemplateNet} < \textbf{DacqGroupChannelEnumTemplateNet}, \textbf{Dacq} \leftarrow \textbf{GroupChannelEnumTemplate}, \textbf{CDeviceGroupChannelInfoTemplateNet} >$	121
CDigOutStimulatorFunctionNet	128
CFilterConfigurationNet	137
CFilterConfigurationRegisterNet	139
CGrapheneFunctionNet	171

CIntanMea_FunctionNet	189
CInterfaceboardFunctionNet	193
CInterfaceboard2FunctionNet	191
CMEA2100x256FunctionNet	339
CMcsBusNet	243
CMcsBus_AxisParametersNet	204
CMcsBus_ExtensionNet	206
CMcsBus_FYIExtensionNet	207
CMcsBus_MotorControlNet	209
CMcsBus_SensorNet	226
CMcsBus_TempSensorNet	236
CMcsBus_VoltageModeNet	238
CMcsUsbDeviceStatePushFunctionNet	297
CMultiwellCallbackFunctionNet	382
CSCUFunctionNet	502
CMeFunctionNet	373
CMeaAudioFunctionNet	341
CMeaDigitalDataFunctionNet	361
CMeaFeedbackFunctionNet	363
CMultiwellOptoStimFunctionNet	391
CPPCFunctionNet	424
CPPS_FunctionNet	434
CProgramPressureCurveNet	439
CPulseGeneratorFunctionNet	440
CRFFunctionNet	448
CRobo_FYIProgram_FunctionNet	453
CRobo_FYITemp_FunctionNet	454
CStimulusFunctionNet	573
CTEERFunctionNet	599
CUsbDeviceConfigurationFunctionNet	611
CW2100_StimulatorFunctionNet	622
CWarnerUssingFunctionNet	630

CWarnerValveControllerDeviceTesterFunctionNet	671
CWirelessBaseFunctionNet	679
CW2100_FunctionNet	614
CWClassicFunctionNet	674
CMcsUsbFunctionPointerContainer	309
CMcsUsbListEntryNet	309
CMcsUsbListNet	315
CMcsUsbNet	318
CExternDTesterDeviceNet	133
CFluidControlDeviceNet	143
CGenericDevelopDeviceNet	150
CGilsonDeviceNet	168
CMcsUsbDacqNet	247
CMeaDeviceNet	353
CMeaUSBDeviceNet	372
CCMOSMeaDeviceNet	113
CGrapheneASICDeviceNet	170
CHLADacqNet	184
CLIH3DeviceNet	195
CMultiwellDeviceNet	384
CWarnerUssingDeviceNet	629
COctoPotDeviceNet	397
CRoboDacqNet	456
CMcsUsbDeviceStatePushNet	298
CWarnerValveControllerDeviceNet	646
CMcsUsbFactoryNet	299
CMeaCleanDeviceNet	344
CMeaCoatDeviceNet	348
CMealmpedanceDeviceNet	368
CMeaSwitchDeviceNet	370
CChannelTestDeviceNet	100
CMultiBatteryChargerDeviceNet	374

CNF_GenDeviceNet	396
COkuvisionStimulatorDeviceNet	401
CPPCDeviceNet	423
CPPS_DeviceNet	433
CPathIdentDeviceNet	405
CPedoterDeviceNet	406
CPeristalticPumpDeviceNet	408
CPgaDeviceNet	409
CPositionIIDeviceNet	411
CPositionImpDeviceNet	420
CRadioControledDevicesNet	443
CRetinaLedDeviceNet	446
CRoboDeviceNet	471
CEncapsulatorDeviceNet	132
CFYIDeviceNet	149
CHLADeviceNet	185
CHiClampDeviceNet	183
CMeasureTableDeviceNet	369
CPPSDeviceNet	438
CPatchServerDeviceNet	404
CRoboStatorDeviceNet	492
CRobolnjectDeviceNet	490
CRoboocyte2DeviceNet	491
CTEERMachineDeviceNet	610
CRoboFluidDeviceNet	487
CSafeISDeviceNet	498
CSerialPortNet	519
CStg200xBasicNet	520
CStg200xDownloadBasicNet	559
CStg200xDownloadNet	566
CSw2to64DeviceNet	584
CTcxDeviceNet	586

CMcsUsbPointerContainer	339
CCMOSMeaDeviceNet::CRegionOfInterestRect	445
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	613
DeviceIdNet	680
DigitalSource< digitalsourceenum >	682
DigitalSourceGeneral	683
DriverVersionNet Exception	684
CUsbExceptionNet	612
FirmwareDestinationNames	690
HeadstageIDTypeObject	696
HeadStageIDTypeState IComparable	698
HeadStageIDType	693
mkfilterNet	699
CRoboDeviceNet::RoboMainLowLevelCommands	702
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	709
CFilterCoefficientsNet::s_FilterAttributesNet	709
CMeaAudioFunctionNet::s_setaudionet	710
CStimulusFunctionNet::SidebandData	711
StgStatusNet stgstreaming	712
CStg200xBasicNet	520
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	712
CMcsUsbFactoryNet::TargetSelectionInfoNet	714
usbSetupPacket_t	714
W2100_StimulusParametersNet	715

9 Class Index

9.1 Class List

Here are the classes, structs, unions and interfaces with brief descriptions:

CW2100_FunctionNet::AudioChannelsNet 96

9.1 Class List

BatteryState	96
BesselFilterBandPassNet	97
BesselFilterHighPassNet	97
BesselFilterLowPassNet	98
ButterworthFilterBandPassNet	99
ButterworthFilterHighPassNet	99
ButterworthFilterLowPassNet	100
CChannelTestDeviceNet CChannelTestDeviceNet	100
CCMOSMea_FunctionNet	102
CCMOSMeaDeviceNet	113
CCreateFilterNet CCreateFilterNet CCreateFilterNet	116
CDacCalibrationFunctionNet 118	
CDacqGroupChannelGenericSelectionNet	120
CDacqGroupChannelSelectionNet	121
CDacqGroupChannelSelectionTemplateNet< DacqGroupChannelEnumTemplateNet, DacqGroupCh 121	annelEnumTemplate
CDeviceGroupChannelInfoGenericNet	124
CDeviceGroupChannelInfoMEA2100_256Net	124
CDeviceGroupChannelInfoNet	125
CDeviceGroupChannelInfoSCUNet	126
CDeviceGroupChannelInfoTemplateNet< DacqGroupChannelEnumTemplateNet >	126
CDeviceGroupChannelInfoW2100Net	127
CDigOutStimulatorFunctionNet CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class	128
CEncapsulatorDeviceNet CEncapsulatorDeviceNet is the to control the MCS HiClamp device	132
CExternDTesterDeviceNet CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)	133
CFilterCoefficientsNet CFILTER CONTROL	135
CFilterConfigurationNet CFILTER CONTROL CONTRO	137
CFilterConfigurationRegisterNet CFILTER CONTROL CONTRO	139
CFilterPropertyNet CFilterProper	141
CFluidControlDeviceNet CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device	143

CFYIDeviceNet CFYIDeviceNet is the class to control the MCS FYI device	149
	143
CGenericDevelopDeviceNet CGenericDevelopDeviceNet is the class to use during development of a new device	150
CGilsonDeviceNet	
CGilsonDeviceNet is the class to control a Gilson device	168
CGrapheneASICDeviceNet	170
CGrapheneFunctionNet CGrapheneFunctionNet is the class to control Graphene device functions	171
CHiClampDeviceNet CHiClampDeviceNet is the to control the MCS HiClamp device	183
CHLADacqNet	184
CHLADeviceNet CHLADeviceNet is the to control the MCS HLA device	185
CMcsUsbDacqNet::CHWInfo Class to provide hardware information about the device	186
CIntanMea_FunctionNet	189
CInterfaceboard2FunctionNet CInterfaceboard2FunctionNet is the class to control the Interfaceboard	191
CInterfaceboardFunctionNet CInterfaceboardFunctionNet is the class to control the Interfaceboard	193
CLIH3DeviceNet CLIH3DeviceNet is the class to access the HEKA LIH3 device	195
CMcsBus_AxisParametersNet	204
CMcsBus_ExtensionNet	206
CMcsBus_FYIExtensionNet	207
CMcsBus_MotorControlNet	209
CMcsBus_SensorNet	226
CMcsBus_TempSensorNet	236
CMcsBus_VoltageModeNet	238
CMcsBusNet	243
CMcsUsbDacqNet Base class for data acquisition devices	247
CMcsUsbDeviceStatePushFunctionNet	297
CMcsUsbDeviceStatePushNet	298
CMcsUsbFactoryNet	299
CMcsUsbFunctionNet	308

9.1 Class List

CMcsUsbFunctionPointerContainer	309
CMcsUsbListEntryNet McsUsbListEntryNet identifies a connected device	309
CMcsUsbListNet Class to handle a list of connected MCS USB devices	315
CMcsUsbNet Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class	318
CMcsUsbPointerContainer	339
CMEA2100_256DacqGroupChannelSelectionNet	339
CMEA2100x256FunctionNet CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include " Stg200xNet.h" to resolve documentation reference	339
CMeaAudioFunctionNet	341
CMeaCleanDeviceNet CMeaCleanDeviceNet is the class to access the MEA Clean device	344
CMeaCoatDeviceNet CMeaCoatDeviceNet is the class to access the MEA Coat device	348
CMeaDeviceNet Base class for MEA data acquisition devices	353
CMeaDigitalDataFunctionNet	361
CMeaFeedbackFunctionNet	363
CMealmpedanceDeviceNet	368
CMeasureTableDeviceNet CMeasureTableDeviceNet is the to control the MCS HLA device	369
CMeaSwitchDeviceNet The class to control the USB-MEA-Switch	370
CMeaUSBDeviceNet Class for data acquisition via ME and MEA USB amplifiers	372
CMeFunctionNet 373	
CMultiBatteryChargerDeviceNet CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device	374
CMultiwellCallbackFunctionNet CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator	382
CMultiwellDeviceNet CMultiwellDeviceNet is the class to access the Multiwell device	384
CMultiwellOptoStimFunctionNet CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device	391

CNF_GenDeviceNet	396
COctoPotDeviceNet	397
COkuvisionStimulatorDeviceNet	401
CPatchServerDeviceNet CPatchServerDeviceNet is the class to control the MCS PatchServer device	404
CPathIdentDeviceNet	405
CPedoterDeviceNet 406	
CPeristalticPumpDeviceNet CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump	408
CPgaDeviceNet	409
CPositionIIDeviceNet CPositionIIDeviceNet is the class to control PositionII devices	411
CPositionImpDeviceNet CPositionImpDeviceNet is the class to access the Position/Imp devices	420
CPPCDeviceNet	423
CPPCFunctionNet CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump	424
CPPS_DeviceNet	433
CPPS_FunctionNet	434
CPPSDeviceNet CPPS4plus1DeviceNet is the to control the MCS HLA device	438
CProgramPressureCurveNet CProgramPressureCurveNet is the class to program pressure curves	439
CPulseGeneratorFunctionNet CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking	440
CRadioControledDevicesNet	443
CCMOSMeaDeviceNet::CRegionOfInterestRect	445
CRetinaLedDeviceNet	446
CRFFunctionNet CRFFunctionNet is the class to control RF devices	448
CRobo_FYIProgram_FunctionNet	453
CRobo_FYITemp_FunctionNet	454
CRoboDacqNet	456
CRoboDeviceNet CRoboDeviceNet is the base class for all Robo platform based devices	471
CRoboFluidDeviceNet	487

9.1 Class List 21

CRobolnjectDeviceNet CRobolnjectDeviceNet is the to control the MCS Robolnject device	490
CRoboocyte2DeviceNet CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device	491
CRoboStatorDeviceNet	492
CSafeISDeviceNet 498	
CSCUDacqGroupChannelSelectionNet	501
CSCUFunctionNet CSCUFunctionNet is the class to control the SCU device	502
CSerialPortNet CSerialPortNet	519
CStg200xBasicNet Base class for the Stg200x	520
CStg200xDownloadBasicNet CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device	559
CStg200xDownloadNet Main class for the STG download mode This class implements the STG download mode interface.	566
CStimulusFunctionNet	573
CSw2to64DeviceNet The class to control the MCS-USB-Sw2to64 device	584
CTcxDeviceNet Class to control a Temperature Controller (TCX)	586
CTEERFunctionNet CTEERFunctionNet is the class to control the TEER device	599
CTEERMachineDeviceNet	610
CUsbDeviceConfigurationFunctionNet CUsbDeviceConfigurationFunctionNet is the class to configure the USB firmware	611
CUsbExceptionNet Exception class that is thrown in case of an USB error	612
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	613
CW2100_FunctionNet	614
CW2100_StimulatorFunctionNet	622
CW2100DacqGroupChannelSelectionNet	628
CWarnerUssingDeviceNet CWarnerUssingDeviceNet is the class to control the Ussing device	629
CWarnerUssingFunctionNet CWarnerUssingFunctionNet is the class to control the Ussing device	630

CWarnerValveControllerDeviceNet CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller	646
CWarnerValveControllerDeviceTesterFunctionNet	
CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the	
Warner Valve Controller Device Tester	671
CWClassicFunctionNet	674
CWirelessBaseFunctionNet	679
DeviceIdNet	
Device Id	680
DigitalSource< digitalsourceenum >	682
DigitalSourceGeneral	683
DriverVersionNet	
Class gives firmware versions of the device's firmware destinations	684
FirmwareDestinationNames	690
HeadStageIDType	693
HeadstageIDTypeObject	696
HeadStageIDTypeState	698
mkfilterNet	699
CRoboDeviceNet::RoboMainLowLevelCommands	702
CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands	709
CFilterCoefficientsNet::s_FilterAttributesNet	709
CMeaAudioFunctionNet::s_setaudionet	710
CStimulusFunctionNet::SidebandData	711
StgStatusNet	712
CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData	712
CMcsUsbFactoryNet::TargetSelectionInfoNet	714
usbSetupPacket_t	714
W2100_StimulusParametersNet	715

10 Namespace Documentation

10.1 Mcs Namespace Reference

Namespaces

• Usb

10.2 Mcs::Usb Namespace Reference

Classes

- class CChannelTestDeviceNet
- class CCMOSMea_FunctionNet
- class CCMOSMeaDeviceNet
- class CCreateFilterNet
- class ButterworthFilterLowPassNet
- · class ButterworthFilterHighPassNet
- · class ButterworthFilterBandPassNet
- class BesselFilterLowPassNet
- · class BesselFilterHighPassNet
- · class BesselFilterBandPassNet
- class CDeviceGroupChannelInfoTemplateNet
- class CDeviceGroupChannelInfoGenericNet
- class CDeviceGroupChannelInfoNet
- class CDeviceGroupChannelInfoW2100Net
- class CDeviceGroupChannelInfoSCUNet
- class CDeviceGroupChannelInfoMEA2100 256Net
- class CDacqGroupChannelSelectionTemplateNet
- class CDacqGroupChannelGenericSelectionNet
- class CDacqGroupChannelSelectionNet
- class CW2100DacqGroupChannelSelectionNet
- · class CSCUDacqGroupChannelSelectionNet
- class CMEA2100 256DacgGroupChannelSelectionNet
- · class CDacCalibrationFunctionNet
- class CDigOutStimulatorFunctionNet

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

class CExternDTesterDeviceNet

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

· class CGrapheneFunctionNet

CGrapheneFunctionNet is the class to control Graphene device functions

· class CInterfaceboard2FunctionNet

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

class CInterfaceboardFunctionNet

CInterfaceboardFunctionNet is the class to control the Interfaceboard

· class CLIH3DeviceNet

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

class CMEA2100x256FunctionNet

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

class CMeaCleanDeviceNet

CMeaCleanDeviceNet is the class to access the MEA Clean device.

class CMeaCoatDeviceNet

CMeaCoatDeviceNet is the class to access the MEA Coat device.

· class CMultiBatteryChargerDeviceNet

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

· class CMultiwellCallbackFunctionNet

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

· class CMultiwellDeviceNet

CMultiwellDeviceNet is the class to access the Multiwell device.

class CMultiwellOptoStimFunctionNet

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

- class CPedoterDeviceNet
- · class CPositionIIDeviceNet

CPositionIIDeviceNet is the class to control PositionII devices

class CPositionImpDeviceNet

CPositionImpDeviceNet is the class to access the Position/Imp devices

class CPPCFunctionNet

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

class CPulseGeneratorFunctionNet

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

class CRFFunctionNet

CRFFunctionNet is the class to control RF devices

· class CSCUFunctionNet

CSCUFunctionNet is the class to control the SCU device

class CTEERFunctionNet

CTEERFunctionNet is the class to control the TEER device

class CUsbDeviceConfigurationFunctionNet

CUsbDeviceConfigurationFunctionNet is the class to configure the USB firmware

class CWarnerUssingDeviceNet

CWarnerUssingDeviceNet is the class to control the Ussing device

· class CWarnerUssingFunctionNet

CWarnerUssingFunctionNet is the class to control the Ussing device

class CWarnerValveControllerDeviceNet

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

class CWarnerValveControllerDeviceTesterFunctionNet

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

struct DeviceIdNet

Device Id.

- · class CFilterCoefficientsNet
- · class CFilterConfigurationNet
- · class CFilterConfigurationRegisterNet
- class CFilterPropertyNet
- class CFluidControlDeviceNet

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

· class CGenericDevelopDeviceNet

CGenericDevelopDeviceNet is the class to use during development of a new device.

• class CGilsonDeviceNet

CGilsonDeviceNet is the class to control a Gilson device.

- class CGrapheneASICDeviceNet
- · class CIntanMea FunctionNet
- class CMcsBusNet
- · class CMcsBus_MotorControlNet
- class CMcsBus VoltageModeNet
- · class CMcsBus_AxisParametersNet
- · class CMcsBus_SensorNet
- class CMcsBus_TempSensorNet
- · class CMcsBus_ExtensionNet
- class CMcsBus_FYIExtensionNet
- class CSerialPortNet
- class usbSetupPacket t
- class CMcsUsbDeviceStatePushFunctionNet

- class CMcsUsbDeviceStatePushNet
- · class CMcsUsbFactoryNet
- · class CMcsUsbFunctionPointerContainer
- class CMcsUsbFunctionNet
- · class CMcsUsbListEntryNet

McsUsbListEntryNet identifies a connected device.

class CMcsUsbListNet

Class to handle a list of connected MCS USB devices.

class CUsbExceptionNet

Exception class that is thrown in case of an USB error.

- class FirmwareDestinationNames
- · class DriverVersionNet

Class gives firmware versions of the device's firmware destinations.

- class CMcsUsbPointerContainer
- class CMcsUsbNet

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

- · class DigitalSourceGeneral
- · class DigitalSource
- class StgStatusNet
- · class CMeaAudioFunctionNet
- · class CMeaDeviceNet

Base class for MEA data acquisition devices.

class CMeaUSBDeviceNet

Class for data acquisition via ME and MEA USB amplifiers

- · class CMeaDigitalDataFunctionNet
- class CMeaFeedbackFunctionNet
- class CMealmpedanceDeviceNet
- · class CMeaSwitchDeviceNet

The class to control the USB-MEA-Switch.

- class CMeFunctionNet
- class mkfilterNet
- · class CNF GenDeviceNet
- · class COctoPotDeviceNet
- · class COkuvisionStimulatorDeviceNet
- · class CPathIdentDeviceNet
- class CPeristalticPumpDeviceNet

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

- class CPgaDeviceNet
- · class CPPCDeviceNet
- · class CPPS DeviceNet
- class CPPS FunctionNet
- · class CProgramPressureCurveNet

CProgramPressureCurveNet is the class to program pressure curves

- class CRadioControledDevicesNet
- · class CRetinaLedDeviceNet
- class CRobo_FYITemp_FunctionNet
- class CRobo_FYIProgram_FunctionNet
- class CRoboDacqNet
- · class CHLADacqNet
- class CRoboDeviceNet

CRoboDeviceNet is the base class for all Robo platform based devices

· class CRoboStatorDeviceNet

· class CRoboocyte2DeviceNet

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

class CRobolnjectDeviceNet

CRobolnjectDeviceNet is the to control the MCS Robolnject device

· class CHiClampDeviceNet

CHiClampDeviceNet is the to control the MCS HiClamp device

· class CEncapsulatorDeviceNet

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

class CHLADeviceNet

CHLADeviceNet is the to control the MCS HLA device

class CPPSDeviceNet

CPPS4plus1DeviceNet is the to control the MCS HLA device

· class CMeasureTableDeviceNet

CMeasureTableDeviceNet is the to control the MCS HLA device

class CFYIDeviceNet

CFYIDeviceNet is the class to control the MCS FYI device

class CPatchServerDeviceNet

CPatchServerDeviceNet is the class to control the MCS PatchServer device

- class CTEERMachineDeviceNet
- class CRoboFluidDeviceNet
- · class CSafeISDeviceNet
- · class CStg200xDownloadNet

Main class for the STG download mode This class implements the STG download mode interface.

class CStg200xBasicNet

Base class for the Stg200x.

class CStg200xDownloadBasicNet

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

- class CStimulusFunctionNet
- · class CSw2to64DeviceNet

The class to control the MCS-USB-Sw2to64 device.

· class CTcxDeviceNet

Class to control a Temperature Controller (TCX)

class CMcsUsbDacqNet

Base class for data acquisition devices.

- struct W2100 StimulusParametersNet
- class HeadStageIDType
- class HeadStageIDTypeState
- · class HeadstageIDTypeObject
- class BatteryState
- class CW2100_StimulatorFunctionNet
- class CW2100 FunctionNet
- class CWClassicFunctionNet
- class CWirelessBaseFunctionNet

Enumerations

```
    enum class enCMosMeaChipType {
        unknown = 0 ,
        nMos16LV = 1 ,
        nMos32LV = 3 ,
        nMos36LN = 6 ,
        nMos64LN = 7 }
```

```
    enum class DeviceEnumNet {

 MCS DEVICE ANY,
 MCS_GENERIC_DEVELOPMENT DEVICE,
 MCS_DEVICE_USB,
 MCS_MCCARD_DEVICE,
 MCS STG DEVICE,
 MCS MC STIMULUS DEVICE,
 MCS MEAUSB DEVICE,
 MCS MEA DEVICE,
 MCS OCTOPOT DEVICE,
 MCS_TERSENS_DEVICE,
 MCS_PGA_DEVICE,
 MCS_PCX_DEVICE,
 MCS TCX DEVICE,
 MCS_FCX_DEVICE,
 MCS_RETINA_LED_DEVICE,
 MCS MEA SWITCH DEVICE.
 MCS MEA IMPEDANCE DEVICE,
 MCS_CHANNELTEST_DEVICE,
 MCS_SW2TO64_DEVICE,
 MCS RETINA AMS DONGLE,
 MCS PATHIDENT DEVICE,
 MCS_ROBO_DEVICE,
 MCS_ROBOOCYTE2_DEVICE,
 MCS ROBOINJECT DEVICE,
 MCS_HICLAMP_DEVICE,
 MCS_PATCHSERVER_DEVICE,
 MCS ENCAPSULATOR DEVICE,
 MCS MEASURETABLE DEVICE,
 MCS FYI DEVICE,
 MCS_HLA_DEVICE
 MCS_PPS_DEVICE,
 MCS_PPS5_DEVICE,
 MCS_OKUVISION_STIMULATOR_DEVICE,
 MCS_NF_GEN_DEVICE,
 MCS_SAFEIS_DEVICE,
 MCS_PERISTALTIC_PUMP_DEVICE,
 MCS_EXTERN_BC_TESTER_DEVICE,
 MCS_EXTERN_D_TESTER_DEVICE,
 MCS SOFTWARE DONGLE DEVICE,
 MCS_MEA_CLEAN_DEVICE,
 MCS_MEA_COAT_DEVICE,
 MCS_SMARTIMPLANT_DEVICE,
 MCS MBC08 DEVICE,
 MCS PEDOTER DEVICE,
 MCS_PPC_DEVICE,
 WARNER_VALVE_CONTROL_DEVICE = 7000,
 WARNER USSING DEVICE,
 HEKA LIH3 DEVICE = 8000,
 ALA_VC3_DEVICE = 9990,
 MCS_DEVICE_USB_CYPRESS = 9991 }
    Enumerates the group of MCS devices to connect to.

    enum class VendorldEnumNet {

 Any = -1,
 None = 0,
 MCS = MCS VENDOR ID,
 PCI = 0x10E8
 Cypress = CYPRESS VENDOR ID,
```

```
ALA_VC3 = ALA_VC3_VENDOR_ID }
    Enumerates the group of MCS devices to connect to.

    enum class ProductIdEnumNet {

 Any = -1,
 None = 0.
 LegacyMeaUsb = MCS_PRODUCT_ID_MEAUSB ,
 ALA VC3 = ALA VC3 VENDOR ID,
 Cypress FX1 = CY FX1 PRODUCT ID,
 Cypress FX2 = CY FX2 PRODUCT ID,
 Cypress_FX3 = CY_FX3_PRODUCT_ID,
 MC Card = MCS PRODUCT ID MC CARD,
 Campden Ci4600EphysVideoDataIntegrator = MCS PRODUCT ID CAMPDEN CI4600EPHYS VIDEO↔
  DATA INTEGRATOR,
 HekaLIH30 = MCS PRODUCT ID HEKA LIH30,
 HekaEPC10Single = MCS_PRODUCT_ID_HEKA_EPC10_SINGLE,
 HekaEPC10Double = MCS PRODUCT ID HEKA EPC10 DOUBLE,
 HekaEPC10Triple = MCS PRODUCT ID HEKA EPC10 TRIPLE,
 HekaEPC10Quadro = MCS_PRODUCT_ID_HEKA_EPC10_QUADRO,
 HekaLIH406 = MCS_PRODUCT_ID_HEKA_LIH_406,
 HekaLIH816 = MCS_PRODUCT_ID_HEKA_LIH_816,
 HekaITEV100 = MCS_PRODUCT_ID_HEKA_ITEV_100,
 HekaPG610 = MCS_PRODUCT_ID_HEKA_PG_610,
 HekaPG611 = MCS_PRODUCT_ID_HEKA_PG_611,
 HekaPG612 = MCS PRODUCT ID HEKA PG 612,
 HekaPG618 = MCS_PRODUCT_ID_HEKA_PG_618,
 HekaPG690 = MCS PRODUCT ID HEKA PG 690,
 HekaEPCLite = MCS PRODUCT ID HEKA EPC Lite,
 STG = MCS PRODUCT ID STG,
 Octopot = MCS_PRODUCT_ID_OCTOPOT,
 Tersens = MCS PRODUCT ID TERSENS,
 Dotriapot = MCS PRODUCT ID DOTRIAPOT,
 HLA = MCS PRODUCT ID HLA,
 STG400x = MCS_PRODUCT_ID_STG400x ,
 STG4002 = MCS PRODUCT ID STG4002,
 STG4004 = MCS PRODUCT ID STG4004,
 STG4008 = MCS PRODUCT ID STG4008,
 STG400x_opto = MCS_PRODUCT_ID_STG400x_OPTO ,
 STG4002 opto = MCS PRODUCT ID STG4002 OPTO,
 STG4004 opto = MCS PRODUCT ID STG4004 OPTO,
 STG4008_opto = MCS_PRODUCT_ID_STG4008 OPTO,
 STG5 = MCS_PRODUCT_ID_STG5 ,
 STG3008_FA = MCS_PRODUCT_ID_STG3008 FA,
 MultiwellOptoStim = MCS_PRODUCT_ID_MULTIWELLOPTOSTIM,
 Generic = MCS_PRODUCT_ID_GENERIC,
 PGA = MCS PRODUCT ID PGA,
 PCX = MCS PRODUCT ID PCX,
 TCX = MCS PRODUCT ID TCX.
 FCX = MCS PRODUCT ID FCX,
 FCB = MCS PRODUCT ID FCB,
 TC01 = MCS PRODUCT ID TC01,
 TC02 = MCS PRODUCT ID TC02,
 Retina_LED = MCS_PRODUCT_ID_RETINA_LED ,
 AMS Dongle = MCS PRODUCT ID RETINA AMS DONGLE,
 Okuvision Stimulator = MCS PRODUCT ID OKUVISION STIMULATOR,
 ExternBCTester = MCS PRODUCT ID RETINAIMPLANT EXTERNBCTESTER,
 Triggerbox_IMS = MCS_PRODUCT_ID_RIAG_TRIGGERBOX_IMS,
 Triggerbox AMS = MCS PRODUCT ID RIAG TRIGGERBOX AMS
 Triggerbox AMS3 = MCS PRODUCT ID RIAG TRIGGERBOX AMS3,
```

```
ExternDTester = MCS_PRODUCT_ID_RETINAIMPLANT_EXTERNDTESTER,
FunkDongleS = MCS PRODUCT ID RIAG FUNKDONGLES
ExternSTester = MCS_PRODUCT_ID_RIAG_EXTERNSTESTER ,
DongleS = MCS_PRODUCT_ID_RIAG_DONGLES,
Triggerbox R5 = MCS PRODUCT ID RIAG TRIGGERBOX R5,
MEA Switch = MCS PRODUCT ID MEA SWITCH,
MEA Impedance = MCS PRODUCT ID MEA IMPEDANCE,
ChannelTest = MCS PRODUCT ID CHANNELTEST,
Sw2to64 = MCS PRODUCT ID SW2TO64,
PeristalticPump = MCS_PRODUCT_ID_PERISTALTIC_PUMP,
MEA_Switch_2_1 = MCS_PRODUCT_ID_MEA_SWITCH_2_1,
MEA_Switch_4_2 = MCS_PRODUCT_ID_MEA_SWITCH_4_2,
PPS4plus1 = MCS_PRODUCT_ID_PPS4plus1 ,
PPS5 = MCS PRODUCT ID PPS5,
PPS2 = MCS_PRODUCT_ID_PPS2,
PPS5_DIG = MCS_PRODUCT_ID_PPS5_DIG ,
MEA Clean = MCS PRODUCT ID MEA CLEAN,
MEA Coat = MCS PRODUCT ID MEA COAT,
Multiwell_ICC = MCS_PRODUCT_ID_MULTIWELL_ICC ,
MBC08 = MCS PRODUCT ID MBC08,
PPC = MCS PRODUCT ID PPC,
MEA1060 = MCS PRODUCT ID MEA1060,
MEA_Sanofi = MCS_PRODUCT_ID_MEA_SANOFI,
ME256 = MCS PRODUCT ID ME256,
ME128 = MCS PRODUCT ID ME128,
ME64 = MCS_PRODUCT_ID_ME64,
ME32 = MCS_PRODUCT_ID_ME32,
ME16 = MCS PRODUCT ID ME16,
MEA2100 Mini Usb develop = MCS PRODUCT ID MEA2100 MINI USB DEVELOP,
MEA256 = MCS PRODUCT ID MEA256,
MEA2100 = MCS_PRODUCT_ID_MEA2100
MEA2100_32 = MCS_PRODUCT_ID_MEA2100_32,
MEA2100 Lite = MCS PRODUCT ID MEA21 LITE,
Multiwell = MCS_PRODUCT_ID_MULTIWELL,
MEA2100_256 = MCS_PRODUCT_ID_MEA2100_256,
ME2100 = MCS_PRODUCT_ID_ME2100,
MEA2100BetaScreen = MCS PRODUCT ID MEA2100 BETA SCREEN,
MEA2100_Mini = MCS_PRODUCT_ID_MEA2100_MINI,
TBSI Dacq = MCS PRODUCT ID TBSI DACQ,
Multiwell MEA Mini = MCS PRODUCT ID MULTIWELL MEA MINI,
Whole Cell Patch = MCS PRODUCT ID WHOLE CELL PATCH,
eCube = MCS PRODUCT ID ECUBE,
Graphene ASIC = MCS PRODUCT ID GRAPHENE ASIC,
GE2100 = MCS PRODUCT ID GE2100,
Multiboot = MCS PRODUCT ID MULTIBOOT,
WPA8 = MCS_PRODUCT_ID_WPA8,
WPA4 = MCS_PRODUCT_ID_WPA4,
WPA16 = MCS PRODUCT ID WPA16,
WPA32 = MCS PRODUCT ID WPA32,
W2100 = MCS_PRODUCT_ID_W2100,
NeuroChip = MCS PRODUCT ID NEUROCHIP,
UsbTest = MCS PRODUCT ID USB TEST,
SoftwareDongle = MCS PRODUCT ID SOFTWAREDONGLE,
PathIdent = MCS_PRODUCT_ID_PATHIDENT,
NF_Gen = MCS_PRODUCT_ID_NF_GEN ,
SafelS = MCS PRODUCT ID SAFEIS,
Encapsulator = MCS_PRODUCT_ID_ENCAPSULATOR,
NeurochipConfig = MCS_PRODUCT_ID_NEUROCHIP_CONFIG,
```

```
MeasureTable = MCS_PRODUCT_ID_MEASURETABLE,
 Robooycte2 = MCS PRODUCT ID ROBOOCYTE2,
 Robolnject = MCS_PRODUCT_ID_ROBOINJECT,
 HiClamp = MCS_PRODUCT_ID_HICLAMP ,
 PatchServer = MCS PRODUCT ID PATCHSERVER,
 Dilutor = MCS PRODUCT ID DILUTOR,
 HiClamp4Uart = MCS PRODUCT ID HICLAMP4UART,
 IM16S16KRA = MCS PRODUCT ID IM16S16KRA,
 IM64KRB = MCS PRODUCT ID IM64KRB,
 IS32KRA = MCS PRODUCT ID IS32KRA,
 IM64KRC = MCS PRODUCT ID IM64KRC
 IM16S8KRA = MCS_PRODUCT_ID_IM16S8KRA ,
 IM16KRC = MCS_PRODUCT_ID_IM16KRC ,
 SmartImplant = MCS PRODUCT ID SMARTIMPLANT,
 PositionImp = MCS_PRODUCT_ID_POSITION_IMP,
 PositionBase = MCS_PRODUCT_ID_POSITION_BASE
 PositionIICentralUnit = MCS_PRODUCT_ID_POSITIONII_CENTRAL_UNIT,
 PositionIIBase = MCS PRODUCT ID_POSITIONII_BASE,
 GrapheneProjectTestDevice = MCS PRODUCT ID GRAPHENE PROJECT TEST DEVICE,
 Pos900 = MCS PRODUCT ID POS900,
 Neptun = MCS PRODUCT ID NEPTUN,
 Warner Valve Control = MCS PRODUCT ID WARNER VALVE CONTROL,
 Warner_TEER_Machine = MCS_PRODUCT_ID_WARNER_TEER_MACHINE,
 Warner_Ussing = MCS_PRODUCT_ID_WARNER_USSING }
    Enumerates the group of MCS devices to connect to.

    enum class McsBusTypeEnumNet {

 MCS ANY BUS = -1,
 MCS UNDEFINED BUS = 0,
 MCS USB BUS,
 MCS PCI BUS }
    Enumerates the bus to use, either USB, PCI or any

    enum class McsUsbSpeedEnumNet {

 LowSpeed = 0,
 FullSpeed = 1,
 HighSpeed = 2
 SuperSpeed = 3,
 UnknownSpeed = 0xff }
    Enumerates the current connection speed of the device

    enum class CFirmwareDestinationNet {

 FPGA NORMAL = 0,
 DSP = MCSUSB DEST DSP,
 USB = MCSUSB DEST USB,
 MCU1 = MCSUSB DEST MCU1 .
 MCSBUS1 = MCSUSB DEST MCSBUS1,
 MCSBUS2 = MCSUSB_DEST_MCSBUS2,
 MCSBUS3 = MCSUSB_DEST_MCSBUS3,
 MCSBUS4 = MCSUSB_DEST_MCSBUS4,
 MCSBUS5 = MCSUSB DEST MCSBUS5,
 MCSBUS6 = MCSUSB_DEST_MCSBUS6,
 MCSBUS7 = MCSUSB_DEST_MCSBUS7,
 MCSBUS8 = MCSUSB DEST MCSBUS8,
 MCSBUS9 = MCSUSB DEST MCSBUS9.
 MCSBUS10 = MCSUSB DEST MCSBUS10,
 MCSBUS11 = MCSUSB DEST MCSBUS11,
 MCSBUS12 = MCSUSB DEST MCSBUS12,
 MCSBUS13 = MCSUSB DEST MCSBUS13.
 MCSBUS14 = MCSUSB_DEST_MCSBUS14,
 MCSBUS15 = MCSUSB_DEST_MCSBUS15,
```

```
MCSBUS0 = MCSUSB_DEST_MCSBUS0,
BUSNUMBERO = MCSUSB DEST BUSNUMBERO,
BUSOMCSBUS1 = MCSUSB_DEST_BUSO_MCSBUS1,
BUSOMCSBUS2 = MCSUSB_DEST_BUS0_MCSBUS2,
BUSOMCSBUS3 = MCSUSB DEST BUSO MCSBUS3,
BUSOMCSBUS4 = MCSUSB DEST BUSO MCSBUS4,
BUSOMCSBUS5 = MCSUSB DEST BUSO MCSBUS5,
BUSOMCSBUS6 = MCSUSB DEST BUSO MCSBUS6,
BUSOMCSBUS7 = MCSUSB DEST BUSO MCSBUS7,
BUSOMCSBUS8 = MCSUSB DEST BUSO MCSBUS8,
BUS0MCSBUS9 = MCSUSB_DEST_BUS0_MCSBUS9,
BUSOMCSBUS10 = MCSUSB_DEST_BUS0_MCSBUS10,
BUSOMCSBUS11 = MCSUSB_DEST_BUSO_MCSBUS11,
BUSOMCSBUS12 = MCSUSB DEST BUSO MCSBUS12,
BUSOMCSBUS13 = MCSUSB_DEST_BUS0_MCSBUS13,
BUSOMCSBUS14 = MCSUSB_DEST_BUS0_MCSBUS14,
BUSOMCSBUS15 = MCSUSB DEST BUSO MCSBUS15,
BUSOMCSBUSO = MCSUSB DEST BUSO MCSBUSO,
BUSNUMBER1 = MCSUSB_DEST_BUSNUMBER1,
BUS1MCSBUS1 = MCSUSB DEST BUS1 MCSBUS1,
BUS1MCSBUS2 = MCSUSB DEST BUS1 MCSBUS2,
BUS1MCSBUS3 = MCSUSB DEST BUS1 MCSBUS3,
BUS1MCSBUS4 = MCSUSB_DEST_BUS1_MCSBUS4,
BUS1MCSBUS5 = MCSUSB_DEST_BUS1_MCSBUS5,
BUS1MCSBUS6 = MCSUSB DEST BUS1 MCSBUS6,
BUS1MCSBUS7 = MCSUSB_DEST_BUS1_MCSBUS7,
BUS1MCSBUS8 = MCSUSB_DEST_BUS1_MCSBUS8,
BUS1MCSBUS9 = MCSUSB DEST BUS1 MCSBUS9,
BUS1MCSBUS10 = MCSUSB DEST BUS1 MCSBUS10,
BUS1MCSBUS11 = MCSUSB DEST BUS1 MCSBUS11,
BUS1MCSBUS12 = MCSUSB_DEST_BUS1_MCSBUS12,
BUS1MCSBUS13 = MCSUSB_DEST_BUS1_MCSBUS13,
BUS1MCSBUS14 = MCSUSB DEST BUS1 MCSBUS14,
BUS1MCSBUS15 = MCSUSB_DEST_BUS1_MCSBUS15,
BUS1MCSBUS0 = MCSUSB_DEST_BUS1_MCSBUS0,
BUSNUMBER2 = MCSUSB_DEST_BUSNUMBER2,
BUS2MCSBUS1 = MCSUSB DEST BUS2 MCSBUS1,
BUS2MCSBUS2 = MCSUSB_DEST_BUS2_MCSBUS2,
BUS2MCSBUS3 = MCSUSB DEST BUS2 MCSBUS3,
BUS2MCSBUS4 = MCSUSB DEST BUS2 MCSBUS4,
BUS2MCSBUS5 = MCSUSB DEST BUS2 MCSBUS5,
BUS2MCSBUS6 = MCSUSB_DEST_BUS2_MCSBUS6,
BUS2MCSBUS7 = MCSUSB DEST BUS2 MCSBUS7,
BUS2MCSBUS8 = MCSUSB DEST BUS2 MCSBUS8,
BUS2MCSBUS9 = MCSUSB DEST BUS2 MCSBUS9,
BUS2MCSBUS10 = MCSUSB_DEST_BUS2_MCSBUS10,
BUS2MCSBUS11 = MCSUSB_DEST_BUS2_MCSBUS11,
BUS2MCSBUS12 = MCSUSB DEST BUS2 MCSBUS12,
BUS2MCSBUS13 = MCSUSB DEST BUS2 MCSBUS13,
BUS2MCSBUS14 = MCSUSB_DEST_BUS2_MCSBUS14,
BUS2MCSBUS15 = MCSUSB DEST BUS2 MCSBUS15,
BUS2MCSBUS0 = MCSUSB DEST BUS2 MCSBUS0,
PIC = MCSUSB DEST PIC,
PIC2 = MCSUSB_DEST_PIC2,
PIC3 = MCSUSB_DEST_PIC3,
PIC4 = MCSUSB_DEST_PIC4,
PIC5 = MCSUSB_DEST_PIC5,
PIC6 = MCSUSB_DEST_PIC6,
```

```
PIC7 = MCSUSB DEST PIC7,
PIC8 = MCSUSB DEST PIC8,
PIC9 = MCSUSB_DEST_PIC9,
PIC10 = MCSUSB_DEST_PIC10,
PIC11 = MCSUSB DEST PIC11,
PIC12 = MCSUSB DEST PIC12,
ChannelPIC = MCSUSB DEST CHANNELPIC.
Bootstrap = MCSUSB DEST BOOTSTRAP,
BootstrapOtherCypress = MCSUSB DEST BOOTSTAP OTHER CYPRESS,
ALTERA = MCSUSB DEST ALTERA,
FPGA2 = MCSUSB DEST FPGA2,
FPGA3 = MCSUSB_DEST_FPGA3,
FPGA4 = MCSUSB_DEST_FPGA4,
FPGA5 = MCSUSB DEST FPGA5,
FPGA6 = MCSUSB_DEST_FPGA6,
FPGA7 = MCSUSB_DEST_FPGA7,
FPGA8 = MCSUSB DEST FPGA8,
FPGA9 = MCSUSB DEST FPGA9.
FPGA10 = MCSUSB DEST FPGA10,
FPGA11 = MCSUSB DEST FPGA11,
FPGA12 = MCSUSB DEST FPGA12,
FPGA13 = MCSUSB DEST FPGA13,
FPGA14 = MCSUSB_DEST_FPGA14,
FPGA15 = MCSUSB DEST FPGA15,
FPGA16 = MCSUSB DEST FPGA16,
FPGA GOLD = XILINX DEST GOLDEN,
ALTERA_GOLD = (MCSUSB_DEST_ALTERA | XILINX_DEST_GOLDEN),
FPGA2 GOLD = (MCSUSB DEST FPGA2 | XILINX DEST GOLDEN),
FPGA3 GOLD = (MCSUSB DEST FPGA3 | XILINX DEST GOLDEN),
FPGA4 GOLD = (MCSUSB DEST FPGA4 | XILINX DEST GOLDEN),
FPGA5_GOLD = (MCSUSB_DEST_FPGA5 | XILINX_DEST_GOLDEN),
FPGA6_GOLD = (MCSUSB_DEST_FPGA6 | XILINX_DEST_GOLDEN),
FPGA7 GOLD = (MCSUSB DEST FPGA7 | XILINX DEST GOLDEN),
FPGA8_GOLD = (MCSUSB_DEST_FPGA8 | XILINX_DEST_GOLDEN),
FPGA9_GOLD = (MCSUSB_DEST_FPGA9 | XILINX_DEST_GOLDEN),
FPGA10_GOLD = (MCSUSB_DEST_FPGA10 | XILINX_DEST_GOLDEN),
FPGA11 GOLD = (MCSUSB DEST FPGA11 | XILINX DEST GOLDEN),
FPGA12_GOLD = (MCSUSB_DEST_FPGA12 | XILINX_DEST_GOLDEN),
FPGA13 GOLD = (MCSUSB DEST FPGA13 | XILINX DEST GOLDEN),
FPGA14 GOLD = (MCSUSB DEST FPGA14 | XILINX DEST GOLDEN),
FPGA15 GOLD = (MCSUSB DEST FPGA15 | XILINX DEST GOLDEN),
FPGA16 GOLD = (MCSUSB DEST FPGA16 | XILINX DEST GOLDEN),
FPGA BASE = XILINX DEST BASEIMAGE,
ALTERA BASE = (MCSUSB DEST ALTERA | XILINX DEST BASEIMAGE),
FPGA2 BASE = (MCSUSB DEST FPGA2 | XILINX DEST BASEIMAGE).
FPGA3_BASE = (MCSUSB_DEST_FPGA3 | XILINX_DEST_BASEIMAGE),
FPGA4 BASE = (MCSUSB DEST FPGA4 | XILINX DEST BASEIMAGE),
FPGA5 BASE = (MCSUSB DEST FPGA5 | XILINX DEST BASEIMAGE),
FPGA6 BASE = (MCSUSB DEST FPGA6 | XILINX DEST BASEIMAGE),
FPGA7_BASE = (MCSUSB_DEST_FPGA7 | XILINX_DEST_BASEIMAGE),
FPGA8 BASE = (MCSUSB DEST FPGA8 | XILINX DEST BASEIMAGE),
FPGA9 BASE = (MCSUSB DEST FPGA9 | XILINX DEST BASEIMAGE),
FPGA10 BASE = (MCSUSB DEST FPGA10 | XILINX DEST BASEIMAGE),
FPGA11_BASE = (MCSUSB_DEST_FPGA11 | XILINX_DEST_BASEIMAGE),
FPGA12_BASE = (MCSUSB_DEST_FPGA12 | XILINX_DEST_BASEIMAGE),
FPGA13 BASE = (MCSUSB DEST FPGA13 | XILINX DEST BASEIMAGE),
FPGA14_BASE = (MCSUSB_DEST_FPGA14 | XILINX_DEST_BASEIMAGE),
FPGA15_BASE = (MCSUSB_DEST_FPGA15 | XILINX_DEST_BASEIMAGE),
```

```
FPGA16_BASE = (MCSUSB_DEST_FPGA16 | XILINX_DEST_BASEIMAGE),
 FPGA BOOTSTRAP = XILINX DEST BOOTSTRAP
 ALTERA BOOTSTRAP = (MCSUSB DEST ALTERA | XILINX DEST BOOTSTRAP),
 DEST_TARGET1 = FLASH_DEST_TARGET1 ,
 DEST_TARGET2 = FLASH_DEST_TARGET2 ,
 DEST TARGET3 = FLASH DEST TARGET3,
 DEST TARGET4 = FLASH DEST TARGET4.
 DEST TARGET5 = FLASH DEST TARGET5,
 DEST TARGET6 = FLASH DEST TARGET6,
 DEST TARGET7 = FLASH DEST TARGET7,
 DEST_TARGET8 = FLASH_DEST_TARGET8 ,
 DEST_TARGET9 = FLASH_DEST_TARGET9 ,
 DEST_TARGET10 = FLASH_DEST_TARGET10 ,
 DEST TARGET11 = FLASH DEST TARGET11,
 DEST_TARGET12 = FLASH_DEST_TARGET12,
 DEST_TARGET13 = FLASH_DEST_TARGET13,
 DEST TARGET14 = FLASH DEST TARGET14,
 DEST_TARGET15 = FLASH_DEST_TARGET15,
 DEST_TARGET_MASK = FPGA_DEST_TARGET_MASK ,
 DEST FX3 TARGET MASK = FX3 DEST TARGET MASK.
 ALTERA TARGET1 = (MCSUSB DEST ALTERA | FLASH DEST TARGET1),
 ALTERA TARGET2 = (MCSUSB DEST ALTERA | FLASH DEST TARGET2),
 ALTERA_TARGET3 = (MCSUSB_DEST_ALTERA | FLASH_DEST_TARGET3),
 USB TARGET1 = (MCSUSB DEST USB | FLASH DEST TARGET1),
 USB TARGET2 = (MCSUSB DEST USB | FLASH DEST TARGET2),
 USB TARGET3 = (MCSUSB DEST USB | FLASH DEST TARGET3),
 UnknownDest = MCSUSB_DEST_UNKNOWN }
    Enumerates the destination processor for the firmware.

    enum class DigitalTargetEnumNet {

 Digout = (MEA COMMAND << 16) + MEA MEA21 DIGOUT SOURCE,
 Digstream = (MEA COMMAND << 16) + MEA MEA21 DIGSTREAM SOURCE,
 DacqTrigger = (MEA COMMAND << 16) + MEA MEA21 DACQTRIGGER SOURCE,
 StgTrigger = (STG200x_COMMAND << 16) + STG200x_TRIGGER_SOURCE,
 StgListModeTrigger = (STG200x COMMAND << 16) + STG200x MEA21 LISTMODE TRIGGERSOURCE
 DigOutStimulatorStartTrigger = (MEA COMMAND << 16) + MEA DIGOUT STG START TRIGGER ←
 SOURCE,
 DigOutStimulatorStopTrigger = (MEA_COMMAND << 16) + MEA_DIGOUT_STG_STOP_TRIGGER_←
 SOURCE,
 DigStreamToReceiver = (MEA COMMAND << 16) + MEA DIGSTREAMTORECEIVER SOURCE }
    Enumerates the Digital Targets for Digital Sources

    enum class MEA2100DigitalSourceEnumNet {

 DigitalInOfOutPort = 0,
 Digitalln = 16,
 DigitalPulse = 32,
 Feedback = 64,
 AuxIn = 96,
 Zero = 98,
 One = 99,
 HS1Trigger1Status = 100,
 HS1Trigger2Status = 102,
 HS1Trigger3Status = 104,
 HS1Trigger4Status = 106.
 HS1Trigger5Status = 108,
 HS1Trigger6Status = 110,
 HS1Sideband1 = 112,
 HS1Sideband2 = 128.
 HS1Sideband3 = 144,
```

```
HS1Sideband4 = 160,
 HS1Sideband5 = 176,
 HS1Sideband6 = 192,
 HS2Trigger1Status = 208,
 HS2Trigger2Status = 210,
 HS2Trigger3Status = 212,
 HS2Trigger4Status = 214,
 HS2Trigger5Status = 216,
 HS2Trigger6Status = 218,
 HS2Sideband1 = 220.
 HS2Sideband2 = 236,
 HS2Sideband3 = 252,
 HS2Sideband4 = 268,
 HS2Sideband5 = 284,
 HS2Sideband6 = 300,
 PulseGenerator = 316,
 DigitalOutStimulator = 320,
 DigitalData = 336,
 DeviceRunStatus = 368,
 LastPosition = 372 }
     Enumerates the digital source of the MEA2100 device.
• enum class MultiwellDigitalSourceEnumNet {
 DigitalInOfOutPort = 0,
 DigitalIn = 16,
 DigitalPulse = 32,
 Feedback = 64,
 AuxIn = 96,
 Zero = 98,
 One = 99,
 HS1Trigger1Status = 100,
 HS1Trigger2Status = 102,
 HS1Sideband1 = 112,
 HS1Sideband2 = 128,
 HS2Trigger1Status = 208,
 HS2Trigger2Status = 210,
 HS2Sideband1 = 220,
 HS2Sideband2 = 236,
 PulseGenerator = 316,
 DigitalOutStimulator = 320,
 DigitalData = 336,
 DeviceRunStatus = 368,
 LastPosition = 372 }
     Enumerates the digital source of the Multiwell device.

    enum class CMOSMEA5000DigitalSourceEnumNet {

 DigitalInOfOutPort = 0,
 DigitalIn = 16,
 DigitalPulse = 32,
 Feedback = 64,
 AuxIn = 96,
 Zero = 98,
 One = 99,
 HS1Trigger1Status = 100,
 HS1Trigger2Status = 102,
 HS1Trigger3Status = 104,
 HS1Sideband1 = 112,
 HS1Sideband2 = 128,
 HS1Sideband3 = 144,
 PulseGenerator = 316,
```

```
DigitalOutStimulator = 320,
 DigitalData = 336,
 DeviceRunStatus = 368,
 LastPosition = 372 }
     Enumerates the digital source of the CMOSMEA5000 device.
• enum class W2100DigitalSourceEnumNet {
 DigitalInOfOutPort = 0,
 DigitalIn = 16,
 DigitalPulse = 32,
 Feedback = 64,
 AuxIn = 96,
 Zero = 98,
 One = 99,
 PulseGenerator = 100,
 DigDataFromReceiver = 128,
 DigitalOutStimulator = 192,
 DigitalData = 208,
 DeviceRunStatus = 240,
 DigStreamFromReceiver = 256,
 LastPosition = 320 }
    Enumerates the digital source of the W2100 device.

    enum class SCUDigitalSourceEnumNet {

 DigitalInOfOutPort = (0x00 << 8) + 0,
 DigitalIn = (0x00 << 8) + 16,
 DigitalPulse = (0x01 << 8),
 Feedback = (0x02 << 8),
 AuxIn = (0x03 << 8) + 0,
 Zero = (0x03 << 8) + 2,
 One = (0x03 << 8) + 3,
 PulseGenerator = (0x03 << 8) + 8,
 DigitalOutStimulator = (0x03 << 8) + 16,
 DigitalData = (0x04 << 8),
 DeviceRunStatus = (0x05 << 8) + 0,
 SCU1HS1Sideband1TriggerStatus = (0x40 << 8) + 0.
 SCU1HS1Sideband2TriggerStatus = (0x40 << 8) + 2,
 SCU1HS2Sideband1TriggerStatus = (0x40 << 8) + 4,
 SCU1HS2Sideband2TriggerStatus = (0x40 << 8) + 6
 SCU1HS3Sideband1TriggerStatus = (0x40 << 8) + 8,
 SCU1HS3Sideband2TriggerStatus = (0x40 << 8) + 10,
 SCU1HS4Sideband1TriggerStatus = (0x40 << 8) + 12,
 SCU1HS4Sideband2TriggerStatus = (0x40 << 8) + 14,
 SCU1LEDStim1TriggerStatus = (0x40 << 8) + 16.
 SCU1LEDStim2TriggerStatus = (0x40 << 8) + 18,
 SCU1LEDStim3TriggerStatus = (0x40 << 8) + 20,
 SCU1LEDStim4TriggerStatus = (0x40 << 8) + 22,
 SCU1HS1Sideband1 = (0x42 << 8),
 SCU1HS1Sideband2 = (0x43 << 8),
 SCU1HS2Sideband1 = (0x44 << 8),
 SCU1HS2Sideband2 = (0x45 << 8),
 SCU1HS3Sideband1 = (0x46 << 8),
 SCU1HS3Sideband2 = (0x47 << 8),
 SCU1HS4Sideband1 = (0x48 << 8)
 SCU1HS4Sideband2 = (0x49 << 8).
 SCU1LEDStim1 = (0x4A << 8),
 SCU1LEDStim2 = (0x4B << 8),
 SCU1LEDStim3 = (0x4C << 8),
 SCU1LEDStim4 = (0x4D << 8),
 SCU2HS1Sideband1TriggerStatus = (0x80 << 8) + 0,
```

```
SCU2HS1Sideband2TriggerStatus = (0x80 << 8) + 2,
SCU2HS2Sideband1TriggerStatus = (0x80 << 8) + 4,
SCU2HS2Sideband2TriggerStatus = (0x80 << 8) + 6,
SCU2HS3Sideband1TriggerStatus = (0x80 << 8) + 8,
SCU2HS3Sideband2TriggerStatus = (0x80 << 8) + 10,
SCU2HS4Sideband1TriggerStatus = (0x80 << 8) + 12,
SCU2HS4Sideband2TriggerStatus = (0x80 << 8) + 14,
SCU2LEDStim1TriggerStatus = (0x80 << 8) + 16,
SCU2LEDStim2TriggerStatus = (0x80 << 8) + 18,
SCU2LEDStim3TriggerStatus = (0x80 << 8) + 20.
SCU2LEDStim4TriggerStatus = (0x80 << 8) + 22
SCU2HS1Sideband1 = (0x82 << 8),
SCU2HS1Sideband2 = (0x83 << 8),
SCU2HS2Sideband1 = (0x84 << 8),
SCU2HS2Sideband2 = (0x85 << 8),
SCU2HS3Sideband1 = (0x86 << 8),
SCU2HS3Sideband2 = (0x87 << 8),
SCU2HS4Sideband1 = (0x88 << 8).
SCU2HS4Sideband2 = (0x89 << 8),
SCU2LEDStim1 = (0x8A << 8)
SCU2LEDStim2 = (0x8B << 8)
SCU2LEDStim3 = (0x8C << 8).
SCU2LEDStim4 = (0x8D << 8),
LastPosition = (0xFF << 8)}
   Enumerates the digital source of the SCU device.
enum class MEA2100 256DigitalSourceEnumNet {
DigitalInOfOutPort = (0x00 << 8) + 0,
DigitalIn = (0x00 << 8) + 16,
DigitalPulse = (0x01 << 8),
Feedback = (0x02 << 8),
AuxIn = (0x03 << 8) + 0,
Zero = (0x03 << 8) + 2
One = (0x03 << 8) + 3,
DeviceRunStatus = (0x03 << 8) + 4,
PulseGenerator = (0x03 << 8) + 8.
DigitalOutStimulator = (0x03 << 8) + 16,
DigitalData = (0x04 << 8),
HS1Trigger1Status = (0x40 << 8) + 0,
HS1Trigger2Status = (0x40 << 8) + 2
HS1Trigger3Status = (0x40 << 8) + 4
HS1Trigger4Status = (0x40 << 8) + 6,
HS1Trigger5Status = (0x40 << 8) + 8
HS1Trigger6Status = (0x40 << 8) + 10
HS1Trigger7Status = (0x40 << 8) + 12,
HS1Trigger8Status = (0x40 << 8) + 14
HS1Trigger9Status = (0x40 << 8) + 16
HS1Trigger10Status = (0x40 << 8) + 18.
HS1Trigger11Status = (0x40 << 8) + 20,
HS1Trigger12Status = (0x40 << 8) + 22
HS1Trigger13Status = (0x40 << 8) + 24
HS1Trigger14Status = (0x40 << 8) + 26
HS1Trigger15Status = (0x40 << 8) + 28,
HS1Trigger16Status = (0x40 << 8) + 30,
HS1Trigger17Status = (0x41 << 8) + 0,
HS1Trigger18Status = (0x41 << 8) + 2
HS1Sideband1 = (0x42 << 8),
HS1Sideband2 = (0x43 << 8),
HS1Sideband3 = (0x44 << 8),
```

```
HS1Sideband4 = (0x45 << 8),
 HS1Sideband5 = (0x46 << 8),
 HS1Sideband6 = (0x47 << 8),
 HS1Sideband7 = (0x48 << 8),
 HS1Sideband8 = (0x49 << 8),
 HS1Sideband9 = (0x4A << 8),
 HS1Sideband10 = (0x4B << 8).
 HS1Sideband11 = (0x4C << 8),
 HS1Sideband12 = (0x4D << 8),
 HS1Sideband13 = (0x4E << 8).
 HS1Sideband14 = (0x4F << 8),
 HS1Sideband15 = (0x50 << 8),
 HS1Sideband16 = (0x51 << 8),
 HS1Sideband17 = (0x52 << 8),
 HS1Sideband18 = (0x53 << 8),
 HS2Trigger1Status = (0x80 << 8) + 0,
 HS2Trigger2Status = (0x80 << 8) + 2
 HS2Trigger3Status = (0x80 << 8) + 4
 HS2Trigger4Status = (0x80 << 8) + 6
 HS2Trigger5Status = (0x80 << 8) + 8,
 HS2Trigger6Status = (0x80 << 8) + 10,
 HS2Trigger7Status = (0x80 << 8) + 12
 HS2Trigger8Status = (0x80 << 8) + 14,
 HS2Trigger9Status = (0x80 << 8) + 16,
 HS2Trigger10Status = (0x80 << 8) + 18,
 HS2Trigger11Status = (0x80 << 8) + 20,
 HS2Trigger12Status = (0x80 << 8) + 22,
 HS2Trigger13Status = (0x80 << 8) + 24,
 HS2Trigger14Status = (0x80 << 8) + 26.
 HS2Trigger15Status = (0x80 << 8) + 28,
 HS2Trigger16Status = (0x80 << 8) + 30,
 HS2Trigger17Status = (0x81 << 8) + 0,
 HS2Trigger18Status = (0x81 << 8) + 2
 HS2Sideband1 = (0x82 << 8),
 HS2Sideband2 = (0x83 << 8),
 HS2Sideband3 = (0x84 << 8),
 HS2Sideband4 = (0x85 << 8),
 HS2Sideband5 = (0x86 << 8),
 HS2Sideband6 = (0x87 << 8),
 HS2Sideband7 = (0x88 << 8),
 HS2Sideband8 = (0x89 << 8),
 HS2Sideband9 = (0x8A << 8),
 HS2Sideband10 = (0x8B << 8),
 HS2Sideband11 = (0x8C << 8),
 HS2Sideband12 = (0x8D << 8),
 HS2Sideband13 = (0x8E << 8),
 HS2Sideband14 = (0x8F << 8),
 HS2Sideband15 = (0x90 << 8),
 HS2Sideband16 = (0x91 << 8),
 HS2Sideband17 = (0x92 << 8),
 HS2Sideband18 = (0x93 << 8),
 LastPosition = (0xFF << 8) }
     Enumerates the digital source of the MEA2100-256 device.

    enum class TBSI DACQDigitalSourceEnumNet {

 DigitalInOfOutPort = (0x00 << 8) + 0,
 DigitalIn = (0x00 << 8) + 16,
 DigitalPulse = (0x01 << 8),
 Feedback = (0x02 << 8),
```

```
AuxIn = (0x03 << 8) + 0,
  Zero = (0x03 << 8) + 2
  One = (0x03 << 8) + 3,
  DeviceRunStatus = (0x03 << 8) + 4,
  PulseGenerator = (0x03 << 8) + 8,
  DigitalOutStimulator = (0x03 << 8) + 16,
  DigitalData = (0x04 << 8),
  HS1DigitalData1 = (0x30 << 8),
  HS2DigitalData1 = (0x70 << 8),
  LastPosition = (0xFF << 8)}
     Enumerates the digital source of the TBSI-DACQ device.
• enum class TriggerSourceEnumNet {
  tsNone = 0,
  tsDigitalIn1 = 1,
  tsDigitalIn2 = 2,
  tsDigitalIn3 = 3,
 tsDigitalIn4 = 4,
 tsDigitalIn5 = 5,
 tsDigitalIn6 = 6,
  tsDigitalIn7 = 7,
 tsDigitalIn8 = 8,
 tsDigitalIn9 = 9,
  tsDigitalIn10 = 10,
  tsDigitalIn11 = 11,
 tsDigitalIn12 = 12,
 tsDigitalIn13 = 13,
  tsDigitalIn14 = 14,
  tsDigitalIn15 = 15,
  tsDigitalIn16 = 16,
  tsDigitalIn17 = 17,
  tsDigitalIn18 = 18,
  tsDigitalIn19 = 19,
 tsDigitalIn20 = 20,
  tsDigitalIn21 = 21,
  tsDigitalIn22 = 22,
  tsDigitalIn23 = 23,
  tsDigitalIn24 = 24,
 tsDigitalIn25 = 25,
  tsDigitalIn26 = 26,
 tsDigitalIn27 = 27,
 tsDigitalIn28 = 28,
 tsDigitalIn29 = 29,
  tsDigitalIn30 = 30,
  tsDigitalIn31 = 31,
  tsDigitalIn32 = 32,
  tsFeedback1 = 33,
  tsFeedback2 = 34.
  tsFeedback3 = 35,
  tsFeedback4 = 36,
  tsFeedback5 = 37,
  tsFeedback6 = 38,
  tsFeedback7 = 39,
  tsFeedback8 = 40,
  tsFeedback9 = 41,
  tsFeedback10 = 42,
  tsFeedback11 = 43,
  tsFeedback12 = 44,
  tsFeedback13 = 45,
```

```
tsFeedback14 = 46,
tsFeedback15 = 47,
tsFeedback16 = 48,
tsFeedback17 = 49,
tsFeedback18 = 50,
tsFeedback19 = 51,
tsFeedback20 = 52,
tsFeedback21 = 53,
tsFeedback22 = 54,
tsFeedback23 = 55.
tsFeedback24 = 56,
tsFeedback25 = 57,
tsFeedback26 = 58,
tsFeedback27 = 59,
tsFeedback28 = 60,
tsFeedback29 = 61,
tsFeedback30 = 62,
tsFeedback31 = 63.
tsFeedback32 = 64,
tsAuxIn1 = 65,
tsAuxIn2 = 66,
tsDigitalPuse0 = 67,
tsDigitalPuse1 = 68,
tsDigitalPuse2 = 69,
tsDigitalPuse3 = 70,
tsDigitalPuse4 = 71,
tsDigitalPuse5 = 72,
tsDigitalPuse6 = 73,
tsDigitalPuse7 = 74,
tsDigitalPuse8 = 75,
tsDigitalPuse9 = 76,
tsDigitalPuse10 = 77,
tsDigitalPuse11 = 78,
tsDigitalPuse12 = 79,
tsDigitalPuse13 = 80,
tsDigitalPuse14 = 81,
tsDigitalPuse15 = 82,
tsDigitalPuse16 = 83,
tsDigitalPuse17 = 84,
tsDigitalPuse18 = 85,
tsDigitalPuse19 = 86,
tsDigitalPuse20 = 87,
tsDigitalPuse21 = 88,
tsDigitalPuse22 = 89,
tsDigitalPuse23 = 90,
tsDigitalPuse24 = 91,
tsDigitalPuse25 = 92,
tsDigitalPuse26 = 93,
tsDigitalPuse27 = 94,
tsDigitalPuse28 = 95,
tsDigitalPuse29 = 96,
tsDigitalPuse30 = 97,
tsDigitalPuse31 = 98,
tsTriggered = 99,
tsSidebandBit8 = 100,
tsDACQCy1Dev1Runs = 101,
tsDACQCy1Dev2Runs = 102,
tsDACQCy2Dev1Runs = 103,
```

```
tsDACQCy2Dev2Runs = 104 }
    Enumerates the trigger source of the MEA2100 device.

    enum class AnalogSourceEnumNet {

 AnalogSource HS1,
 AnalogSource HS2.
 AnalogSource IF }
    Enumerates the analog source of the MEA2100 device.
enum class Stg200xTriggerStatusEnumNet {
 Idle = 0,
 Running = 1,
 Finished = 2,
 Armed = 3
    Enumerates the STG download mode trigger status

    enum class Stg3008FilterAmpAmplificationEnumNet {

 Gain1000 = 0,
 Gain500 = 1,
 Gain200 = 2
 Gain100 = 3
    Enumerates the STG3008FA filter amplifier gains

    enum class RetriggerActionEnumNet {

 Stop = STG200x RETRIGGER STOP
 Restart = STG200x_RETRIGGER_RESTART,
 Ignore = STG200x_RETRIGGER_IGNORE,
 Gate = STG200x RETRIGGER GATEMODE,
 Single = STG200x RETRIGGER SINGLE }
    Enumerates possible retrigger actions for STG200x devices.

    enum class Stg200xSegmentFlagsEnumNet {

 None = 0.
 UpdateTrigger = SEGMENTFLAGS_UPDATETRIGGER ,
 DownloadOnly = SEGMENTFLAGS_DOWNLOADONLY,
 TriggerOnly = SEGMENTFLAGS TRIGGERONLY,
 SyncStart = SEGMENTFLAGS SYNCSTART }
    Enumerates Segmentflag options for STG400x devices.

    enum class Stg200xMultiFileSubmodeEnumNet {

 MultiFileMode = MULTIFILESUBMODE MULTIFILEMODE,
 ExtendedMultiFileMode = MULTIFILESUBMODE_EXTENDEDMULTIFILEMODE }
    Enumerates EnableMultifileMode submodes for STG devices.

    enum class Stg200xDigoutModeEnumNet {

 Monitor = STG200x DIGOUTMODE MONITOR,
 Manual = STG200x DIGOUTMODE MANUAL,
 SYNCOUT1 = STG200x DIGOUTMODE SYNCOUT1,
 SYNCOUT2 = STG200x_DIGOUTMODE_SYNCOUT2,
 SYNCOUT3 = STG200x DIGOUTMODE SYNCOUT3,
 SYNCOUT4 = STG200x DIGOUTMODE SYNCOUT4,
 SYNCOUT5 = STG200x DIGOUTMODE SYNCOUT5,
 SYNCOUT6 = STG200x DIGOUTMODE SYNCOUT6,
 SYNCOUT7 = STG200x_DIGOUTMODE SYNCOUT7,
 SYNCOUT8 = STG200x DIGOUTMODE SYNCOUT8 }
    Enumerates the DigoutMode on STG400x devices.

    enum class DigitalStimulatorTriggerSlopeEnumNet {

 Falling = 0,
 Rising = 1 }
    Enumerates start/stop conditions for DigOut/DigStim trigger. /summary>

    enum class DigitalStimulatorTriggerEventEnumNet {

 Start = 0,
 Stop = 1 }
```

```
Enumerates start/stop event for DigOut/DigStim trigger. /summary>
enum class AdapterTypeEnumNet {
 None = 0,
 MEA60 = 1,
 MEA2x60 = 2,
 MEA120 = 3,
 MEA32 = 4,
 MEA2x32 = 5
 Multiwell96 = 6,
 WirelessTestAdapter = 7,
 MEA252 = 8,
 MEA 2 252 2 = 9,
 MEA_2_252_26Well = 10,
 MEA_2_252_2_9Well = 11,
 MEA_2_252_2 Test = 12,
 TBSI 5 = 13,
 TBSI 15 = 14,
 TBSI_31 = 15,
 TBSI_{63} = 16,
 TBSI 127 = 17,
 TBSI_Reserved = 18,
 Ci4600Intan = 20,
 Unknown = ADAPTER_TYPE_UNKOWN,
 NotApplicable = ADAPTER_TYPE_ENUM_NOT_APPLICABLE }
     Enumerates the adapter type of the MEA2100 device.

    enum class MeaLayoutEnumNet {

 mIUnknown = 0,
 mIMEA60 = 1
     Enumerates the MEA layout of the MEA2100 device.

    enum class DataModeEnumNet {

 Unsigned 16bit = 0,
 Unsigned 24bit = 2,
 Unsigned_32bit = 3,
 Signed_16bit = 8,
 Signed_24bit = 10,
 Signed 32bit = 11 }
     Enumerates the data mode of the device, either 16, 24 or 32 bit, can be signed or unsigned.

    enum class SampleSizeNet {

 SampleSize16Unsigned = 2,
 SampleSize16Signed = 2 + 0x100,
 SampleSize24Unsigned = 3,
 SampleSize24Signed = 3 + 0x100,
 SampleSize32Unsigned = 4,
 SampleSize32Signed = 4 + 0x100,
 SampleSize64Unsigned = 8,
 SampleSize64Signed = 8 + 0x100 }
     Enumerates the data format for ChannelBlock functions.

    enum class SampleDstSizeNet {

 SampleDstSize16 = 2,
 SampleDstSize32 = 4 }
     Enumerates the destination data format for ChannelBlock functions.

    enum class TcxDeviceTypeEnumNet {

 Unknown = 0,
 Regular = 1,
 BMI = 2,
 Nanion = 3,
 Warner = 4 }
```

```
Enumerates the type of TCX devices.
enum class TcxSensorTypeEnumNet {
 Reserved5 = 0,
 Reserved4 = 1,
 Reserved3 = 2,
 Reserved2 = 3,
 Reserved 1 = 4,
 NTC10K = 5,
 PT1000 = 6.
 PT100 = 7
     Enumerates the sensor types for TCX devices

    enum class STG DestinationEnumNet {

 channeldata voltage,
 channeldata current,
 syncoutdata,
 channeldata_positive_voltage,
 channeldata_positive_current,
 rawdata,
 channeldata_current_own_sync,
 channeldata_positive_current_own_sync,
 channeldata_current_own_boost_gnd_sync,
 channeldata positive current own boost gnd sync,
 channeldata current always boost.
 channeldata_current_always_boost_own_sync }
     Enumerates the destination for STG downloads.
• enum class ElectrodeModeEnumNet {
 emAutomatic = 0,
 emManual = 3 }
    Enumerates the mode of each electrode, can be automatic or manual. In automatic mode, the blanking of the
    electrode is controlled by the sideband signal, in manual mode, the stimulation configuration is independant of the
    sideband signal.

    enum class ElectrodeDacMuxEnumNet {

 Ground = 0
 Stg1 = 1
 Stg2 = 2,
 Stg3 = 3
     Enumerates the setting of the Stimulation DAC Multiplexer.

    enum class DacgGroupChannelEnumNet {

 HeadstageElectrodeGroup = 0x00,
 InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP.
 DSPDataGroup = DSPDATACHANNELSGROUP,
 Headstage1NCBathCurrentGroup = 0x30,
 Headstage1NCCol2CurrentGroup = 0x31,
 Headstage1NChipTempGroup = 0x32,
 STG1DACSignalGroup = 0x38,
 LIH30UserADCGroup = 0x50,
 LIH30TestADCGroup = 0x51,
 LIH30ADCModulesGroup = 0x52.
 IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP,
 STG1SidebandsGroup = 0x90,
 STG1TriggerStatusGroup = 0x91,
 DACQ1DigitalGroup = 0xA0,
 AudioTestChannelGroup = AUDIOTESTCHANNELGROUP,
 PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }
     Enumerates the Channel Groups of Datastream

    enum class W2100DacqGroupChannelEnumNet {

 InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP,
```

```
DSPDataGroup = DSPDATACHANNELSGROUP,
WirelessHeadStageAnalogRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0,
WirelessHeadStageStatusRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1,
WirelessHeadStageAnalogRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0,
WirelessHeadStageStatusRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1,
WirelessHeadStageAnalogRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0,
WirelessHeadStageStatusRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1,
WirelessHeadStageAnalogRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0,
WirelessHeadStageStatusRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1,
WirelessHeadStageAnalogRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0.
WirelessHeadStageStatusRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1,
WirelessHeadStageAnalogRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0,
WirelessHeadStageStatusRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1,
WirelessHeadStageAnalogRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0,
WirelessHeadStageStatusRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1,
WirelessHeadStageAnalogRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0,
WirelessHeadStageStatusRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1,
WirelessHeadStageGyroDataRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 + 16,
WirelessHeadStageAccDataRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 16,
WirelessHeadStageGyroDataRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 + 16,
WirelessHeadStageAccDataRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 16,
WirelessHeadStageGyroDataRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 + 16,
WirelessHeadStageAccDataRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 16,
WirelessHeadStageGyroDataRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 + 16,
WirelessHeadStageAccDataRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 16,
WirelessHeadStageGyroDataRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 + 16,
WirelessHeadStageAccDataRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 16,
WirelessHeadStageGyroDataRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0 + 16.
WirelessHeadStageAccDataRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 16.
WirelessHeadStageGyroDataRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0 + 16.
WirelessHeadStageAccDataRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 16,
WirelessHeadStageGyroDataRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0 + 16,
WirelessHeadStageAccDataRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 16,
WirelessHeadStageOptoStimCurrentRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 +
0 + 32,
WirelessHeadStageReservedARE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 32
WirelessHeadStageOptoStimCurrentRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 +
WirelessHeadStageReservedARE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 32
WirelessHeadStageOptoStimCurrentRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 +
WirelessHeadStageReservedARE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 32
WirelessHeadStageOptoStimCurrentRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 +
0 + 32,
WirelessHeadStageReservedARE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 32
WirelessHeadStageOptoStimCurrentRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 +
0 + 32.
WirelessHeadStageReservedARE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 32
WirelessHeadStageOptoStimCurrentRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 +
0 + 32,
WirelessHeadStageReservedARE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 32
WirelessHeadStageOptoStimCurrentRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 +
```

```
0 + 32,
 WirelessHeadStageReservedARE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 32
 WirelessHeadStageOptoStimCurrentRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 +
 0 + 32.
 WirelessHeadStageReservedARE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 32
 WirelessHeadStageReservedBRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 0 + 48
 WirelessHeadStageReservedCRE1HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 0 + 1 + 48
 WirelessHeadStageReservedBRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 0 + 48
 WirelessHeadStageReservedCRE1HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 2 + 1 + 48
 WirelessHeadStageReservedBRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 0 + 48
 WirelessHeadStageReservedCRE1HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 4 + 1 + 48
 WirelessHeadStageReservedBRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 0 + 48
 WirelessHeadStageReservedCRE1HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 0 + 6 + 1 + 48
 WirelessHeadStageReservedBRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 0 + 48
 WirelessHeadStageReservedCRE2HS1 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 0 + 1 + 48
 WirelessHeadStageReservedBRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 0 + 48
 WirelessHeadStageReservedCRE2HS2 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 2 + 1 + 48
 WirelessHeadStageReservedBRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 0 + 48
 WirelessHeadStageReservedCRE2HS3 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 4 + 1 + 48
 WirelessHeadStageReservedBRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 0 + 48
 WirelessHeadStageReservedCRE2HS4 = WIRELESSHEADSTAGEANALOGGROUPBASE + 8 + 6 + 1 + 48
 IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP,
 AudioTestChannelGroup = AUDIOTESTCHANNELGROUP,
 PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }
    Enumerates the W2100 Channel Groups of Datastream

    enum class SCUDacqGroupChannelEnumNet {

 SCU1ElectrodeGroupHS1 = 0x00,
 SCU1ElectrodeGroupHS2 = 0x01,
 SCU1ElectrodeGroupHS3 = 0x02.
 SCU1ElectrodeGroupHS4 = 0x03.
 SCU2ElectrodeGroupHS1 = 0x08,
 SCU2ElectrodeGroupHS2 = 0x09,
 SCU2ElectrodeGroupHS3 = 0x0A.
 SCU2ElectrodeGroupHS4 = 0x0B,
 InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP,
 STG1DACSignalGroup = 0x40,
 STG2DACSignalGroup = 0x41,
 DSPAnalogGroup = DSPDATACHANNELSGROUP,
 DSPDigitalGroup = 0xA0,
 IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP,
```

```
STG1TriggerStatusGroup = 0x90,
 STG1SidebandsGroup = 0x91.
 STG2TriggerStatusGroup = 0x98,
 STG2SidebandsGroup = 0x99,
 AudioTestChannelGroup = AUDIOTESTCHANNELGROUP,
 PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }
     Enumerates the SCU Channel Groups of Datastream

    enum class MEA2100_256DacqGroupChannelEnumNet {

 HS1ElectrodeGroup = 0x00,
 HS2ElectrodeGroup = 0x08,
 InterfaceADCGroup = INTERFACEANALOGCHANNELSGROUP,
 STG1DACSignalGroup = 0x40,
 STG2DACSignalGroup = 0x41,
 DSPAnalogGroup = DSPDATACHANNELSGROUP,
 DSPDigitalGroup = 0xA0.
 IFDigChannelsGroup = INTERFACEDIGITALCHANNELSGROUP,
 STG1TriggerStatusGroup = 0x90,
 STG1SidebandsGroup = 0x91.
 STG2TriggerStatusGroup = 0x98,
 STG2SidebandsGroup = 0x99,
 AudioTestChannelGroup = AUDIOTESTCHANNELGROUP,
 PacketFrameContextGroup = PACKETFRAMECONTEXTGROUP }
    Enumerates the MEA2100-256 Channel Groups of Datastream
 enum class DacqMeaGroupTypeEnumNet {
 AnalogGroup = ANALOG GROUP ,
 DigitalGroup = DIGITAL GROUP.
 FrameContextGroup = FRAME CONTEXT GROUP }
     Enumerations of CMOS MEA Groups to detect wether it is an Analog, Digital or Frame Context Group

    enum class CMOSMeaValueUnitEnumNet {

 NoUnit = 0x00,
 NanoVolt = 0x11.
 PicoAmpere = 0x21
 NanoAmpere = 0x22,
 MicroAmpere = 0x23.
 MilliDegreeCelsius = 0x31 }
     Enumerations of CMOS MEA Units of Values in Data stream

    enum class CMOSMeaInterfaceADCEnumNet {

 IFChannel1 = 0x01,
 IFChannel2 = 0x02,
 IFChannel3 = 0x04,
 IFChannel4 = 0x08,
 IFChannel5 = 0x10,
 IFChannel6 = 0x20,
 IFChannel7 = 0x40,
 IFChannel8 = 0x80 }
     Enumerations of CMOS MEA IF Analog Channels Group Bitmask

    enum class CMOSMeaHeadstage1NCBathCurrentEnumNet { NCBathCurrent = 0x01 }

     Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

    enum class CMOSMeaHeadstage1NCCol2CurrentEnumNet { NCCol2Current = 0x01 }

     Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask

    enum class CMOSMeaHeadstage1NChipTempEnumNet { NChipTemperature = 0x01 }

     Enumerations of CMOS MEA HS Temperature Monitoring Channels Group Bitmask

    enum class CMOSMeaSTG1DACSignalEnumNet {

 DAC1Channel = 0x01,
 DAC2Channel = 0x02.
 DAC3Channel = 0x04,
 DAC4Channel = 0x08 }
```

```
Enumerations of CMOS MEA DAC Stimulation Channels Group Bitmask

    enum class CMOSMealFDigChannelEnumNet {

 DigitalMux = 0x01,
 DigitalInPort = 0x02,
 DigitalOutReg = 0x04,
 FeedbackReg = 0x08,
 DigitalReg = 0x10,
 AuxPort = 0x20 }
     Enumerations of CMOS MEA IF Digital Channels Group Bitmask
• enum class CMOSMeaHS1SidebandEnumNet {
 SBSVector1 = 0x01,
 SBSVector2 = 0x02,
 SBSVector3 = 0x04,
 SBSVector4 = 0x08 }
     Enumerations of CMOS MEA HS STG Sideband Channels Group Bitmask

    enum class CMOSMeaHS1TriggerStatusEnumNet {

  TriggerStatus1 = 0x01,
 TriggerStatus2 = 0x02,
 TriggerStatus3 = 0x04
 TriggerStatus4 = 0x08 }
     Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask
enum class AnalogUnitEnumNet {
 Unknown,
 Volt,
 Ampere,
 Kelvin }

    enum class CMOSMeaPacketFrameContextGroupEnumNet {

 SOFAndCTRLword = 0x01,
 ChecksumAndPacketCounter = 0x02,
 Timestamp = 0x04,
 EOFAndCRC = 0x08 }
     Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

    enum class CMOSMeaBathModeEnumNet {

 Ground = 0x02,
 Stimulation = 0x01,
 CurrentMeasure = 0x00 }
    Enumerations of CMOS MEA Bath Mode
• enum class PatchServAdcModeEnumNet {
 Normal = 0.
 CatchAmp = 1 }

    enum class RoboCurrentModeEnumNet {

 Off = ROBO_CURRENT_OFFMODE,
 Break = ROBO_CURRENT_BREAKMODE,
 Standby = ROBO CURRENT STANDBYMODE,
 Reference = ROBO_CURRENT_REFERENCEMODE,
 Movement = ROBO_CURRENT_MOVEMENTMODE }

    enum class TeerClampModeEnumNet {

 ClampModeVoltage = 0,
 ClampModeCurrent = 1,
 ClampModeOpen = 2,
 ClampModeInternalCalibration = 3 }

    enum class TeerWaveformEnumNet {

 Rectangle = 0,
 Sine = 1 }

    enum class UssingClampModeEnumNet {

 VoltageClamp = 1,
```

```
CurrentClamp = 2,
 OpenClamp = 3,
 Standby = 4,
 ElectrodeOffset = 5 }
 enum class UssingUnitEnumNet {
 Volt = 0,
 Ampere = 1,
 State = 2 }
enum class PlateClampEnumNet {
 Close = 0,
 Open = 1,
 Stop = 2

    enum class PlateClampLockEnumNet {

 Lock = 0,
 Unlock = 1 }

    enum class MultiwellPlateTypeEnumNet {

 Plate Dummy = HS PLATETYPE 0,
 Plate_24W700_100FMA = 1,
 Plate_24W030MGA = 2,
 Plate_{72}W500_{100}PMA = 3,
 Plate_72W500_100FMA = 5,
 Plate_24W700_100FMB = HS_PLATETYPE_6,
 Plate_96W700_100FMA = HS_PLATETYPE_7,
 Plate 96W300 80 1152FMA = HS PLATETYPE 33,
 Plate_96W400_80_1152FMB = HS_PLATETYPE_36,
 Plate_24W300_30_1152GBA = HS_PLATETYPE_40,
 Plate 24W700 100FMC = HS PLATETYPE 44,
 Plate 96W700 100FMB = HS PLATETYPE 48,
 Plate_96W700_100GBC = HS_PLATETYPE_49,
 Plate 96W700 100GBD = HS PLATETYPE 51,
 Plate 24W700 100PBA = HS PLATETYPE 60,
 Plate Dummy 126 = HS PLATETYPE 126,
 Plate_24W300_30GMA = HS_PLATETYPE_193,
 Plate 96W700 100GMA = 194,
 Plate 24W300 30GBA = HS PLATETYPE 195,
 Plate 96W700 100GBA = HS PLATETYPE 224,
 Plate_24W300_30GBB = HS_PLATETYPE_232,
 Plate 96W700_100GBB = HS_PLATETYPE_244,
 No Plate = 255 }
 enum class FpgaldEnumNet {
 DeviceNotConnected = FPGA ID NOT CONNECTED,
 Mea2100Interfaceboard = FPGA ID MEA2100 IF,
 Mea2100Headstage = FPGA_ID_MEA2100_HS ,
 Mea2100STG = FPGA_ID_MEA2100_STG,
 MultiwellHeadstage = FPGA ID HS MULTIWELL,
 MultiwellInterfaceboard = FPGA ID IF MULTIWELL.
 TbsiDacqInterfaceboard = FPGA_ID_TBSI_DACQ_IF,
 TbsiDacqHeadstage = FPGA ID TBSI DACQ HS,
 CmosMeaInterfaceboard = FPGA ID CMOS MEA IF,
 CmosMeaHeadstage = FPGA_ID_CMOS_MEA_HS ,
 Mea2100MultiwellIFB2 = FPGA_ID_MEA2100_MW_IFB2,
 Me2100Interfaceboard = FPGA_ID_ME2100_IFB,
 Me2100InvivoSignalCollectorUnit = FPGA ID ME2100 InvivoSCU,
 Me2100InvitroSignalCollectorUnit = FPGA ID ME2100 InvitroSCU,
 Me2100 32XilinxHeadstage = FPGA ID ME2100 32 XILINX HS,
 Me2100 32PICiCE40Headstage = FPGA ID ME2100 32 PIC ICE40 HS,
 Mea2100 256Interfaceboard = FPGA ID MEA2100 256 IF,
 Mea2100 256Headstage = FPGA ID MEA2100 256 HS,
```

```
W2100Interfaceboard = FPGA ID W2100 IF,
 W2100WirelessReceiver = FPGA ID W2100 REC,
 W2100WirelessReceiverAnalog = FPGA ID W2100 REC ANALOG,
 Mea2100Mini60PICiCE40Headstage = FPGA_ID_MEA2100MINI60_PIC_ICE40_HS,
 Mea2100BetaScreenHeadstage = FPGA ID MEA2100BETASCREEN HS,
 Me2100UPA32Headstage = FPGA ID ME2100UPA32 HS,
 MultiwellMiniHeadstage = FPGA ID MULTIWELL MINI HS.
 Mea2100Mini120Headstage = FPGA ID MEA2100MINI120 HS,
 Mea2100Mini60ECP5Headstage = FPGA ID MEA2100MINI60 ECP5 HS,
 eCubeHeadstage = FPGA ID ECUBE HS.
 Me2100Graphene16 32Headstage = FPGA ID ME2100 GRAPHENE 16 32 HS,
 GrapheneASICHeadstage = FPGA_ID_GRAPHENE_ASIC_HS,
 WholeCellPatchHeadstage = FPGA_ID_WHOLE_CELL_PATCH_HS,
 InterfaceBoard2 = FPGA ID INTERFACEBOARD2,
 W2100IFB2 = FPGA_ID_W2100_IFB2,
 CmosmealFB2 = FPGA_ID_CMOS_MEA_IFB2,
 M4M Dock = FPGA ID M4M DOCK,
 M4M ReUse = FPGA ID M4M REUSE.
 Mea2100LiteHeadstage = FPGA ID MEA2100 LITE HS,
 LIH30Interfaceboard = FPGA ID LIH30 USB IF,
 LIH30ADCCtrl = FPGA ID LIH30 ADC CTRL,
 UssingRail = FPGA ID USSING RAIL,
 UssingChamber = FPGA ID USSING CHAMBER,
 IFB2GoldenInterfaceboard = FPGA ID IFB2 GOLDEN,
 IFB30GoldenInterfaceboard = FPGA ID IFB30 GOLDEN,
 DeviceHasNoHeadstage = FPGA ID HAS NO HS }

    enum class HeadstageIdEnumNet {

 DeviceNotConnected = FPGA ID NOT CONNECTED.
 Mea2100 = FPGA ID MEA2100 HS.
 Multiwell = FPGA ID HS MULTIWELL,
 TbsiDacq = FPGA_ID_TBSI_DACQ_HS,
 CmosMea = FPGA ID CMOS MEA HS,
 InvivoSignalCollectorUnit = FPGA ID ME2100 InvivoSCU,
 InvitroSignalCollectorUnit = FPGA ID ME2100 InvitroSCU,
 Mea2100 256 = FPGA ID MEA2100 256 HS,
 W2100WirelessReceiver = FPGA ID W2100 REC.
 W2100WirelessReceiverAnalog = FPGA_ID_W2100_REC_ANALOG,
 Mea2100 Lite = FPGA ID MEA2100 LITE HS,
 LIH30ADCCtrl = FPGA ID LIH30 ADC CTRL,
 DeviceHasNoHeadstage = FPGA ID HAS NO HS }

    enum class SCU HeadstageldEnumNet {

 DeviceNotConnected = FPGA ID NOT CONNECTED.
 Me2100 32Xilinx = FPGA ID ME2100 32 XILINX HS,
 Me2100_32PICiCE40 = FPGA_ID_ME2100_32_PIC_ICE40_HS,
 Mea2100Mini60PICiCE40 = FPGA_ID_MEA2100MINI60_PIC_ICE40_HS,
 Mea2100BetaScreen = FPGA ID MEA2100BETASCREEN HS,
 Me2100UPA32 = FPGA ID ME2100UPA32 HS,
 MultiwellMini = FPGA ID MULTIWELL MINI HS,
 Mea2100Mini120 = FPGA ID MEA2100MINI120 HS,
 Mea2100Mini60ECP5 = FPGA ID MEA2100MINI60 ECP5 HS,
 eCube = FPGA_ID_ECUBE_HS,
 Me2100Graphene16 32 = FPGA ID ME2100 GRAPHENE 16 32 HS,
 GrapheneASIC = FPGA ID GRAPHENE ASIC HS,
 WholeCellPatch = FPGA ID WHOLE CELL PATCH HS,
 DeviceHasNoHeadstage = FPGA_ID_HAS_NO_HS }

    enum class UsbVendorldEnumNet {

 Unknown = -1,
 None = 0,
```

```
Renesas = 0x1912,
 ASMedia = 0x1b21,
 Intel = 0x8086 }
• enum class FilterCalculationDirectionEnumNet {
 DoubleToInt = 0,
 IntToDouble = 1 }
enum class FilterBandEnumNet {
 Unknown = 0.
 Lowpass = 1,
 Highpass = 2 }
enum class FilterFamilyEnumNet {
 Unknown = 0,
 Bessel = 1,
 Butterworth = 2,
 RC = 3
enum class FilterTypeEnumNet {
 Hardware = 0,
 Software = 1 }
enum class FilterAttributeEnumNet {
 PreCommaB = 0,
 PostCommaB = 1,
 CommaPositionB = 2,
 PreCommaA = 3,
 PostCommaA = 4,
 CommaPositionA = 5 }

    enum class AnalogOut DAC Range EnumNet {

 PlusMinus2Comma5Volts = 0,
 PlusMinus5Volts = 1,
 PlusMinus10Volts = 2 }
enum class PP_Pump_Mode_Type_EnumNet {
 Manual = 1,
 Digital = 2,
 Analog = 3 }
• enum class MbcChargingModeEnumNet {
 StorageCharge = 0,
 FullCharge = 1 }
enum class MbcRatedCapacityEnumNet {
 rc30mAh = 0,
 rc100mAh = 1,
 rc200mAh = 2,
 rc300mAh = 3,
 rcGreater300mAh = 4 }

    enum class MbcChannelStateEnumNet {

 csldleNoBattery = 0,
 csldleChargeFinished = 1,
 csCapacityTestPrecharge = 2,
 csCapacityTestDischarge = 3,
 csRefreshBattery = 4,
 csCharge = 5,
 csDischarge = 6,
 csError = 7 }

    enum class PulseGenerator Mode EnumNet {

 Off = 0.
 AlwaysOn = 1,
 Gated_Low_Active = 2,
 Gated High Active = 3 }
• enum class LIH30_ADC_Channel_EnumNet {
 User\_ADC\_0 = 0,
```

```
User\_ADC\_1 = 1,
 User ADC 2 = 2,
 User\_ADC\_3 = 3,
 User\_ADC\_4 = 4,
 Test\_ADC\_EPC10 = 5,
 ModulA ADC0 = 6,
 ModulA ADC1 = 7,
 ModulA\_ADC2 = 8,
 ModulA ADC3 = 9,
 ModulB\_ADC0 = 10,
 ModulB\_ADC1 = 11,
 ModulB\_ADC2 = 12,
 ModulB\_ADC3 = 13,
 ModulC_ADC0 = 14,
 ModulC_ADC1 = 15,
 ModulC\_ADC2 = 16,
 ModulC ADC3 = 17,
 ModulD ADC0 = 18,
 ModulD_ADC1 = 19,
 ModulD_ADC2 = 20,
 ModulD ADC3 = 21 }
• enum class LIH30_DAC_Channel_EnumNet {
 User_DAC_0 = 0,
 User_DAC_1 = 1,
 User_DAC_2 = 2,
 Test_DAC_EPC10 = 3,
 ModulA_DAC0 = 4,
 ModulA DAC1 = 5,
 ModulB DAC0 = 6.
 ModulB_DAC1 = 7,
 ModulC_DAC0 = 8,
 ModulC DAC1 = 9,
 ModulD DAC0 = 10,
 ModuID_DAC1 = 11 }
enum class LIH30_EPC10_Bus_EnumNet {
 A = 0,
 B = 1
enum class W2100_Accel_Gyro_Select_EnumNet {
 Off = 0,
 GyroOnly = 1,
 AccelOnly = 2,
 Both = 3
    enumerates the accelerometer configuration on the W2100 device
enum class WvcValveModeEnumNet {
 Manual = WVC_VALVE_MODE_MANUAL,
 Digital = WVC_VALVE_MODE_DIGITAL,
 Analog = WVC VALVE MODE ANALOG,
 Table = WVC_VALVE_MODE_TABLE }
    enumerates Wvc valve mode

    enum class WvcDisplayModeEnumNet {

 Work = WVC DISPLAY MODE WORK,
 PC = WVC_DISPLAY_MODE_PC,
 Settings = WVC DISPLAY MODE SETTINGS,
 TouchTest = WVC_DISPLAY_MODE_TOUCH_TEST }
    enumerates Wvc display mode

    enum class PortDirectionEnumNet {

 Output = 0,
 Input = 1 }
```

```
enumerates a port direction

    enum class StimulationLayoutConfigurationEnumNet {

 SingleWell = 1,
 SixWell = 2,
 NineWell = 3 }
     enumerates the layout configuration for the MEA2100-256 device

    enum class ReferenceElectrodeSwitchPositionEnumNet {

 off = 0,
 Ref8 = 1,
 Ref16 = 2
 Ref24 = 3,
 Ref32 = 4
     enumerates the possible positions of the reference electrode switch of the ME2100 device

    enum class ReferenceElectrodeModeEnumNet {

 SubtractionOff = 0,
 SubtractFromAllOther = 1,
 SubtractFromReferenceElectrodeOnly = 2,
 SubtractFromAll = 3 }
     enumerates the electrode subtraction modes

    enum class DigitalDatastreamEnableEnumNet {

 None = 0x0000,
 Mux = 0x0001,
 MuxOtherDevice = 0x0002,
 DigitalInReserverd = 0x0004,
 DigitalIn = 0x0008,
 DigitalOut = 0x0010,
 DigitalOutReserved = 0x0020,
 RegisterLow = 0x0040,
 RegisterHigh = 0x0080,
 FeedbackLow = 0x0100,
 FeedbackHigh = 0x0200,
 Aux = 0x0400,
 PeriodicPulse = 0x0800,
 DigOutStim = 0x1000,
 Hs1Digital = 0x00008000,
 Hs1Trigger = 0x00010000,
 Hs1SidebandLow = 0x00020000,
 Hs1SidebandHigh = 0x00040000,
 Hs2Digital = 0x00800000,
 Hs2Trigger = 0x01000000,
 Hs2SidebandLow = 0x020000000,
 Hs2SidebandHigh = 0x04000000 }
     enumerates the streams available as digital datastream

    enum class loVoltageEnumNet {

 Voltage 3V3 = IFB2 IO VOLTAGE 3V3,
 Voltage_5V0 = IFB2_IO_VOLTAGE_5V0 }
     enumerates the I/O Voltages available on the IFB2
enum class EnSTG200x_STATUS {
 OK,
 NOT CONNECTED,
 DEVICE_NOT_FOUND }
```

Functions

public delegate void OnMcsUsbDeviceState (usbSetupPacket t[^] request)

- private delegate void OnMcsUsbDeviceStateCallback (IntPtr pThis, uint32_t size, IntPtr buffer)
- public delegate void OnUpdateFirmwareStatusChange (String^)
- public delegate void OnUpdateFirmwareProgress (int)
- public delegate void OnDeviceArrivalRemoval (CMcsUsbListEntryNet^ entry)

Delegate to show a device arrival or removal.

- public delegate void OnStgPollStatus (unsigned int status, StgStatusNet[^] stgStatusNet, array< int >[^] index_list)
- public delegate void OnMwPollStatus (unsigned int CurrentTemp, unsigned int PlateState, unsigned int SwitchState)
- public delegate void RoboStatusEventDelegate (array< unsigned char >^ buffer)
- public delegate void OnStg200xDataHandler (uint32 t trigger)
- public delegate void OnStg200xErrorHandler ()
- public delegate void OnChannelData (CMcsUsbDacqNet[^] dacq, int CbHandle, int numFrames)
- public delegate void OnError (String^ msg, int action)

10.2.1 Enumeration Type Documentation

10.2.1.1 AdapterTypeEnumNet enum AdapterTypeEnumNet [strong]

Enumerates the adapter type of the MEA2100 device.

None	
MEA60	
MEA2x60	
MEA120	
MEA32	
MEA2x32	
Multiwell96	
WirelessTestAdapter	
MEA252	
MEA_2_252_2	
MEA_2_252_2_6Well	
MEA_2_252_2_9Well	
MEA_2_252_2_Test	
TBSI_5	
TBSI_15	
TBSI_31	
TBSI_63	
TBSI_127	
TBSI_Reserved	
Ci4600Intan	
Unknown	
NotApplicable	_

10.2.1.2 AnalogOut_DAC_Range_EnumNet enum AnalogOut_DAC_Range_EnumNet [strong]

Enumerator

PlusMinus2Comma5Volts	
PlusMinus5Volts	
PlusMinus10Volts	

10.2.1.3 AnalogSourceEnumNet enum AnalogSourceEnumNet [strong]

Enumerates the analog source of the MEA2100 device.

Enumerator

AnalogSource_HS1	
AnalogSource_HS2	
AnalogSource_IF	

10.2.1.4 AnalogUnitEnumNet enum AnalogUnitEnumNet [strong]

Enumerator

Unknown	
Volt	
Ampere	
Kelvin	

10.2.1.5 **CFirmwareDestinationNet** enum CFirmwareDestinationNet [strong]

Enumerates the destination processor for the firmware.

FPGA_NORMAL	
DSP	The DSP.
USB	The USB controller.
MCU1	The DSP on the MEA2100 system.
MCSBUS1	
MCSBUS2	
MCSBUS3	
MCSBUS4	
MCSBUS5	

MCSBUS6	
MCSBUS7	
MCSBUS8	
MCSBUS9	
MCSBUS10	
MCSBUS11	
MCSBUS12	
MCSBUS13	
MCSBUS14	
MCSBUS15	
MCSBUS0	
BUSNUMBER0	
BUS0MCSBUS1	
BUS0MCSBUS2	
BUS0MCSBUS3	
BUS0MCSBUS4	
BUS0MCSBUS5	
BUS0MCSBUS6	
BUS0MCSBUS7	
BUS0MCSBUS8	
BUS0MCSBUS9	
BUS0MCSBUS10	
BUS0MCSBUS11	
BUS0MCSBUS12	
BUS0MCSBUS13 BUS0MCSBUS14	
BUS0MCSBUS15	
BUS0MCSBUS0	
BUSNUMBER1	
BUS1MCSBUS1	
BUS1MCSBUS2	
BUS1MCSBUS3	
BUS1MCSBUS4	
BUS1MCSBUS5	
BUS1MCSBUS6 BUS1MCSBUS7	
BUS1MCSBUS8	
BUS1MCSBUS9	
BUS1MCSBUS10	
BUS1MCSBUS11 BUS1MCSBUS12	
BUS1MCSBUS13	
BUS1MCSBUS14	
BUS1MCSBUS15	
BUS1MCSBUS0	
BUSNUMBER2	
BUS2MCSBUS1	
BUS2MCSBUS2	
BUS2MCSBUS3	
BUS2MCSBUS4	

DI ION TOODI IO	
BUS2MCSBUS5	
BUS2MCSBUS6	
BUS2MCSBUS7	
BUS2MCSBUS8	
BUS2MCSBUS9	
BUS2MCSBUS10	
BUS2MCSBUS11	
BUS2MCSBUS12	
BUS2MCSBUS13	
BUS2MCSBUS14	
BUS2MCSBUS15	
BUS2MCSBUS0	
PIC	
PIC2	
PIC3	
PIC4	
PIC5	
PIC6	
PIC7	
PIC8	
PIC9	
PIC10	
PIC11	
PIC12	
ChannelPIC	
Bootstrap	
BootstrapOtherCypress	
ALTERA	
FPGA2	
FPGA3	
FPGA3	
FPGA5	
FPGA6	
FPGA7	
FPGA8	
FPGA9	
FPGA10	
FPGA11	
FPGA12	
FPGA13	
FPGA14	
FPGA15	
FPGA16	
FPGA_GOLD	
ALTERA_GOLD	
FPGA2_GOLD	
FPGA3_GOLD	
FPGA4_GOLD	
FPGA5_GOLD	
FPGA6 GOLD	

FPGA7 GOLD	
FPGA8 GOLD	
FPGA9 GOLD	
FPGA10 GOLD	
FPGA11 GOLD	
FPGA12 GOLD	
FPGA13 GOLD	
FPGA13_GOLD	
FPGA14_GOLD	
FPGA15_GOLD	
FPGAT6_GOLD FPGA_BASE	
ALTERA_BASE FPGA2 BASE	
<u> </u>	
FPGA3_BASE	
FPGA4_BASE	
FPGA5_BASE	
FPGA6_BASE	
FPGA7_BASE	
FPGA8_BASE	
FPGA9_BASE	
FPGA10_BASE	
FPGA11_BASE	
FPGA12_BASE	
FPGA13_BASE	
FPGA14_BASE	
FPGA15_BASE	
FPGA16_BASE	
FPGA_BOOTSTRAP	
ALTERA_BOOTSTRAP	
DEST_TARGET1	
DEST_TARGET2	
DEST_TARGET3	
DEST_TARGET4	
DEST_TARGET5	
DEST_TARGET6	
DEST_TARGET7	
DEST_TARGET8	
DEST_TARGET9	
DEST_TARGET10	
DEST_TARGET11	
DEST_TARGET12	
DEST_TARGET14	
DEST_TARGET14	
DEST_TARGET15	
DEST_TARGET_MASK	
DEST_FX3_TARGET_MASK	
ALTERA_TARGET1	
ALTERA_TARGET2	
ALTERA_TARGET3	

1	USB_TARGET1
2	USB_TARGET2
3	USB_TARGET3
st	UnknownDest

10.2.1.6 CMOSMEA5000DigitalSourceEnumNet enum CMOSMEA5000DigitalSourceEnumNet [strong]

Enumerates the digital source of the CMOSMEA5000 device.

Enumerator

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Trigger3Status	
HS1Sideband1	
HS1Sideband2	
HS1Sideband3	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
LastPosition	

10.2.1.7 CMOSMeaBathModeEnumNet enum CMOSMeaBathModeEnumNet [strong]

Enumerations of CMOS MEA Bath Mode

Ground	
Stimulation	
CurrentMeasure	

10.2.1.8 CMOSMeaHeadstage1NCBathCurrentEnumNet enum CMOSMeaHeadstage1NCBathCurrentEnumNet [strong]
Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask
Enumerator
NCBathCurrent
40.04.0. 0110011
10.2.1.9 CMOSMeaHeadstage1NCCol2CurrentEnumNet enum CMOSMeaHeadstage1NCCol2CurrentEnumNet [strong]
Enumerations of CMOS MEA HS Current Monitoring Channels Group Bitmask
Enumerator
NCCol2Current
10.2.1.10 CMOSMeaHeadstage1NChipTempEnumNet enum CMOSMeaHeadstage1NChipTempEnumNet [strong
Enumerations of CMOS MEA HS Temperature Monitoring Channels Group Bitmask
Enumerator
NChipTemperature

10.2.1.11 CMOSMeaHS1SidebandEnumNet enum CMOSMeaHS1SidebandEnumNet [strong]

Enumerations of CMOS MEA HS STG Sideband Channels Group Bitmask

Enumerator

SBSVector1	
SBSVector2	
SBSVector3	
SBSVector4	

10.2.1.12 CMOSMeaHS1TriggerStatusEnumNet enum CMOSMeaHS1TriggerStatusEnumNet [strong]

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

TriggerStatus1	
TriggerStatus2	
TriggerStatus3	
TriggerStatus4	

10.2.1.13 CMOSMealFDigChannelEnumNet enum CMOSMealFDigChannelEnumNet [strong]

Enumerations of CMOS MEA IF Digital Channels Group Bitmask

Enumerator

DigitalMux	
DigitalInPort	
DigitalOutReg	
FeedbackReg	
DigitalReg	
AuxPort	

10.2.1.14 CMOSMeaInterfaceADCEnumNet enum CMOSMeaInterfaceADCEnumNet [strong]

Enumerations of CMOS MEA IF Analog Channels Group Bitmask

Enumerator

IFChannel1	
IFChannel2	
IFChannel3	
IFChannel4	
IFChannel5	
IFChannel6	
IFChannel7	
IFChannel8	

10.2.1.15 CMOSMeaPacketFrameContextGroupEnumNet enum CMOSMeaPacketFrameContextGroupEnumNet [strong]

Enumerations of CMOS MEA HS STG Trigger Status Channels Group Bitmask

SOFAndCTRLword	
ChecksumAndPacketCounter	
Timestamp	
EOFAndCRC	

10.2.1.16 CMOSMeaSTG1DACSignalEnumNet enum CMOSMeaSTG1DACSignalEnumNet [strong]

Enumerations of CMOS MEA DAC Stimulation Channels Group Bitmask

Enumerator

DAC1Channel	
DAC2Channel	
DAC3Channel	
DAC4Channel	

10.2.1.17 CMOSMeaValueUnitEnumNet enum CMOSMeaValueUnitEnumNet [strong]

Enumerations of CMOS MEA Units of Values in Data stream

Enumerator

NoUnit	
NanoVolt	
PicoAmpere	
NanoAmpere	
MicroAmpere	
MilliDegreeCelsius	

10.2.1.18 DacqGroupChannelEnumNet enum DacqGroupChannelEnumNet [strong]

Enumerates the Channel Groups of Datastream

HeadstageElectrodeGroup	
InterfaceADCGroup	
DSPDataGroup	
Headstage1NCBathCurrentGroup	
Headstage1NCCol2CurrentGroup	

Headstage1NChipTempGroup	
STG1DACSignalGroup	
LIH30UserADCGroup	
LIH30TestADCGroup	
LIH30ADCModulesGroup	
IFDigChannelsGroup	
STG1SidebandsGroup	
STG1TriggerStatusGroup	
DACQ1 DigitalGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.19 DacqMeaGroupTypeEnumNet enum DacqMeaGroupTypeEnumNet [strong]

Enumerations of CMOS MEA Groups to detect wether it is an Analog, Digital or Frame Context Group

Enumerator

AnalogGroup	
DigitalGroup	
FrameContextGroup	

10.2.1.20 DataModeEnumNet enum DataModeEnumNet [strong]

Enumerates the data mode of the device, either 16, 24 or 32 bit, can be signed or unsigned.

Enumerator

Unsigned_16bit	
Unsigned_24bit	
Unsigned_32bit	
Signed_16bit	
Signed_24bit	
Signed_32bit	

10.2.1.21 DeviceEnumNet enum DeviceEnumNet [strong]

Enumerates the group of MCS devices to connect to.

MCS_GENERIC_DEVELOPMENT_DEVICE MCS_GENERIC_DEVELOPMENT_DEVICE MCS_MCS_DEVICE_USB MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_MCS_DEVICE MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_		
MCS_MCS_TRIDEUSE MCS_MCS_TRIMULUS_DEVICE MCS_MCS_TRIMULUS_DEVICE MCS_MCS_TRIMULUS_DEVICE Devices which should be accessible from MC_Stimulus. MCS_MCS_MCS_TRIMULUS_DEVICE Devices which should be accessible from MC_Stimulus. MCS_MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_TRIMULUS_DEVICE MCS_MCS_TRIMULUS_DEVICE MCS_MCS_MCS_DEVICE MCS_DEVICE MCS_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_TCX_DEVICE MCS_TCX_DEVICE MCS_TCX_DEVICE MCS_TCX_DEVICE MCS_RETINA_LED_DEVICE MCS_MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_		-
MCS_MCCARD_DEVICE MCS_MCS_STIG_DEVICE MCS_MCS_STIMULUS_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE MCS_MCS_MEAUSB_DEVICE MCS_MCS_MCS_MEAUSC MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_	MCS_GENERIC_DEVELOPMENT_DEVICE	Please use this only for MCS internal development.
MCS_MC_STMULUS_DEVICE MCS_MC_STMULUS_DEVICE MCS_MC_STMULUS_DEVICE MCS_MCS_MEA_DEVICE Connect to an MCS MeaUsb device. MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_MCS_	MCS_DEVICE_USB	To connect to any MCS USB device.
MCS_MC_STIMULUS_DEVICE MCS_MEAUSB_DEVICE MCS_MEAUSB_DEVICE Connect to an MCS MeaUsb device. Connect to an MCS MeaUsb device. MCS_DCTOPOT_DEVICE MCS_TERSENS_DEVICE MCS_PCA_DEVICE MCS_PCA_DEVICE MCS_PCA_DEVICE MCS_PCA_DEVICE Connect to an MCS Tersens device. Connect to an MCS Tersens device. Connect to an MCS Tersens device. Connect to an MCS PCA device. MCS_PCA_DEVICE MCS_TCX_DEVICE MCS_FCA_DEVICE MCS_FCA_DEVICE MCS_METINA_LED_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_CHANNELITEST_DEVICE MCS_CHANNELITEST_DEVICE MCS_RETINA_AMS_DONGLE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE MCS_ROBOO_VEYED_EVICE MCS_ROBOO_VEYED_EVICE MCS_ROBOO_VEYED_EVICE MCS_ROBOON_TEZ_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBOON_TEZ_DEVICE Connect to an MCS Roboon_Verdevice. Connect to an MCS Roboon_Verdevice. MCS_ROBOON_TEZ_DEVICE MCS_ROBOON_TEZ_DEVICE Connect to an MCS Roboon_Verdevice. MCS_ROBOIN_LECT_DEVICE MCS_PATHICLAMP_DEVICE MCS_PATHICLAMP_DEVICE MCS_PATHICLAMP_DEVICE MCS_PATHICLAMP_DEVICE MCS_ROBOON_TEZ_DEVICE Connect to an MCS Roboon_Verdevice. MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE Connect to an MCS Roboon_Verdevice. MCS_PATCHSERVER_DEVICE Connect to an MCS Roboon_Verdevice. MCS_PAS_DEVICE MCS_PAS_DEVICE MCS_PAS_DEVICE Connect to an MCS Roboon_Verdevice. Connect to an MC	MCS_MCCARD_DEVICE	Connect to an MC_Card.
MCS_MEAUSB_DEVICE MCS_MEA_DEVICE Connect to an MCS MeaUsb device. Connect to an MCS MeaUsb device. Connect to an MCS Cotopot device. MCS_PCA_DEVICE MCS_PGA_DEVICE MCS_PGA_DEVICE MCS_PCA_DEVICE MCS_PCA_DEVICE MCS_PCA_DEVICE Connect to an MCS PCX device. Connect to an MCS PCX device. Connect to an MCS PCX device. Connect to an MCS FCX device. MCS_PCX_DEVICE MCS_RETINA_LED_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_SW2T064_DEVICE MCS_SW2T064_DEVICE MCS_RATHIDENT_DEVICE MCS_ROBOINJECT_DEVICE MCS_PATHIDEVICE MCS_PATHIDEVICE	MCS_STG_DEVICE	Connect to an MCS device with STG capability.
MCS_MEA_DEVICE MCS_OCTOPOT_DEVICE MCS_TERSENS_DEVICE MCS_PGA_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_TCX_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_SW2TO64_DEVICE MCS_SW2TO64_DEVICE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_MCS_PATCHSERVER_DEVICE MCS_PTCA_DEVICE MCS_PYLDEVICE MCS_PYLDEVICE MCS_PYLDEVICE MCS_PPS_DEVICE MCS_PRISTALTIC_PUMP_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVI	MCS_MC_STIMULUS_DEVICE	Devices which should be accessible from MC_Stimulus.
MCS_OCTOPOT_DEVICE MCS_TERSENS_DEVICE MCS_PGA_DEVICE MCS_PGA_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_MCS_FCX_DEVICE MCS_MCS_WITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_CHANNELTEST_DEVICE MCS_CHANNELTEST_DEVICE MCS_SATHIDENT_DEVICE MCS_RETINA_AMS_DONGLE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_MCS_PATCHSERVER_DEVICE MCS_PSD_DEVICE MCS_PSD_DEVICE MCS_PSD_DEVICE MCS_PSD_DEVICE MCS_PSD_DEVICE MCS_PSD_DEVICE MCS_PSS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_MCS_MCS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS	MCS_MEAUSB_DEVICE	Connect to an MCS MeaUsb device.
MCS_TERSENS_DEVICE MCS_PCA_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_TCX_DEVICE MCS_TCX_DEVICE MCS_TCX_DEVICE MCS_EX_DEVICE MCS_RETINA_LED_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_MEPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_SW2T064_DEVICE MCS_SW2T064_DEVICE MCS_SW2T064_DEVICE MCS_RETINA_AMS_DONGLE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_RETINA_AMS_DONGLE MCS_RETINA_AMS_DONGLE MCS_ROBO_DEVICE MCS_ROBOO_OPEVICE MCS_ROBOO_OPEVICE MCS_ROBOO_OPEVICE MCS_ROBOO_OPEVICE MCS_ROBOO_OPEVICE MCS_ROBOO_OPEVICE MCS_NOBOO_OPEVICE MCS_HICLAMP_DEVICE MCS_HICLAMP_DEVICE MCS_HICLAMP_DEVICE MCS_MEA_SURTABLE_DEVICE MCS_MEA_SURTABLE_DEVICE MCS_MEA_SURTABLE_DEVICE MCS_MEA_SURTABLE_DEVICE MCS_PYS_DEVICE MCS_PYS_DEVICE MCS_PS_DEVICE MCS_PS_DEVICE MCS_PS_DEVICE MCS_PSS_DEVICE MCS_PSS_DEVICE MCS_SAFEIS_DEVICE MCS_PERISTALTIC_PUMP_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CATA_DEVICE MCS_MEA_CATA_DEVICE MCS_MEA_CATA_DEVICE MCS_MEA_CATA_DEVICE MCS_MEA_CATA_DEVICE MCS_MEA_CATA_DEVICE MCS_SAFEIS_DEVICE MCS_SAFE	MCS_MEA_DEVICE	Connect to an MCS MeaUsb device.
MCS_PGA_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_PCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_FCX_DEVICE MCS_RETINA_LED_DEVICE MCS_RETINA_LED_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_SWITO64_DEVICE MCS_SWITO64_DEVICE MCS_SWITO64_DEVICE MCS_SWITO64_DEVICE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MCS_FTABLE_DEVICE MCS_MCS_FTABLE_DEVICE MCS_PATCHSERVER_DEVICE MCS_MCS_FV_DEVICE MCS_PS_DEVICE MCS_PRIBLIATOR_DEVICE MCS_PS_DEVICE Connect to an MCS PristalticPump device. MCS_PS_CONNECT to an MCS PeristalticPump device. MCS_PS_CONNECT to an MCS PeristalticPump device. MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE Connect to an SeternDTester device. MCS_MEA_CLEAN_DEVICE MCS_PC_DEVICE MCS_PC_DEVICE MCS_PC_DEVICE Connect to a MEA Clean device. MCS_PC_DEVICE MCS_PC_DEVICE MCS_PC_DEVICE Connect to a PCC device. MCS_PC_DEVICE MCS_PC_DEVICE Connect to a PCC device. MCS_PC_DEVICE MCS_PC_DEVICE COnnect to a Warner Valve Control device. MCS_PC_DEVICE	MCS_OCTOPOT_DEVICE	Connect to an MCS Octopot device.
MCS_PCX_DEVICE MCS_TCX_DEVICE Connect to an MCS TCX device. MCS_TCX_DEVICE Connect to an MCS TCX device. MCS_RETINA_LED_DEVICE Connect to an MCS RetineLed device. MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_IMPEDANCE_DEVICE MCS_MEA_IMPEDANCE_DEVICE Connect to an MCS Mea Impedance device. MCS_CHANNELITEST_DEVICE Connect to an MCS Mea Impedance device. MCS_CHANNELITEST_DEVICE Connect to an MCS SW2TO64 device. MCS_RETINA_AMS_DONGLE MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS Dongle (Radio device) MCS_PATHIDENT_DEVICE Connect to an MCS Robo Platform device. MCS_ROBO_DEVICE MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOCYTE2_DEVICE Connect to an MCS Roboocyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolnject device. MCS_HICLAMP_DEVICE Connect to an MCS Robolnject device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MCS_ENCAPSULATOR_DEVICE MCS_HILA_DEVICE Connect to an MCS Encapsulator device. MCS_HILA_DEVICE Connect to an MCS FIV device. MCS_HILA_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS NF-Gen device. MCS_PPS_DEVICE Connect to an MCS PeristalticPump device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternBCTester device. MCS_MEA_CLEAN_DEVICE Connect to an ExternBCTester device. MCS_MEA_CLEAN_DEVICE Connect to an Sextern Device. MCS_MEA_CLEAN_DEVICE Connect to an Sextern Device. MCS_MEA_CLEAN_DEVICE Connect to an MCS PC device. MCS_PEDOTER_DEVICE Connect to a WEAT Clean device. MCS_PEDOTER_DEVICE Connect to a WEAT Clean device. MCS_PEDOTER_DEVICE Connect to a Weather Valve Control device. MCS_PEDOTER_DEVICE Connect to a Warner Valve Control device. WARNE	MCS_TERSENS_DEVICE	Connect to an MCS Tersens device.
MCS_TCX_DEVICE MCS_RETINA_LED_DEVICE Connect to an MCS FCX device. MCS_MEA_SWITCH_DEVICE Connect to an MCS RetineLed device. MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_IMPEDANCE_DEVICE Connect to an MCS Mea Impedance device. MCS_CHANNELTEST_DEVICE MCS_SW2TO64_DEVICE MCS_SW2TO64_DEVICE Connect to an MCS Mea Impedance device. MCS_SW2TO64_DEVICE MCS_SW2TO64_DEVICE MCS_SW2TO64_DEVICE MCS_SW2TO64_DEVICE MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS_Dongle (Radio device) MCS_PATHIDENT_DEVICE Connect to an MCS Retina AMS_Dongle (Radio device) MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOINJECT_DEVICE MCS_ROBOINJECT_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_MEASURETABLE_DEVICE MCS_MCS_PTS_DEVICE MCS_PS_DEVICE MCS_PERISTALTIC_PUMP_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_DETESTER_DEVICE MCS_MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_COAT_DEVICE MCS_PEDOTER_DEVICE MCS_PC_DEVICE MC	MCS_PGA_DEVICE	Connect to an MCS PGA device.
MCS_FCX_DEVICE Connect to an MCS FCX device. MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_IMPEDANCE_DEVICE Connect to an MCS Mea Impedance device. MCS_CHANNELTEST_DEVICE Connect to an MCS Mea Impedance device. MCS_SW2TO64_DEVICE Connect to an MCS SW2TO64 device. MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS Dongle (Radio device) MCS_PATHIDENT_DEVICE Connect to an MCS Pathident device. MCS_ROBO_CVTE2_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOCYTE2_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOLNECT_DEVICE Connect to an MCS Roboolnject device. MCS_HICLAMP_DEVICE Connect to an MCS Roboolnject device. MCS_PATCHSERVER_DEVICE Connect to an MCS PathServer device. MCS_PATCHSERVER_DEVICE Connect to an MCS PathServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_MEASURETABLE_DEVICE Connect to an MCS Encapsulator device. MCS_PYI_DEVICE Connect to an MCS PYI device. MCS_PS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PS_DEVICE Connect to an MCS PPS device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS NF-Gen device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_DC_TESTER_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_DC_TESTER_DEVICE Connect to an MCS PeristalticPump device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_PEDOTER_DEVICE Connect to a PPC device. MCS_PEDOTER_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_C	MCS_PCX_DEVICE	Connect to an MCS PCX device.
MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE MCS_MEA_SWITCH_DEVICE Connect to an MCS Mea Switch device. MCS_MEA_IMPEDANCE_DEVICE MCS_CHANNELTEST_DEVICE Connect to an MCS ChannelTest device. MCS_CHANNELTEST_DEVICE Connect to an MCS ChannelTest device. MCS_SW2T064_DEVICE MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS Dongle (Radio device) MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBO_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOINJECT_DEVICE MCS_ROBOINJECT_DEVICE MCS_HICLAMP_DEVICE Connect to an MCS Robolopiect device. MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_HLA_DEVICE MCS_PS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PS_DEVICE MCS_OKUVISION_STIMULATOR_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_PERISTALTIC_PUMP_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_SAFEIS_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_COAT_DEVICE MCS_MEA_COAT_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SAFEIS_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_COAT_DEVICE MCS_MEA_COAT_DEVICE MCS_PEDOTER_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PEO device. MCS_PEDOTER_DEVICE COnnect to a Warner Valve Control device. MCS_PEDOTER_DEVICE COnnect to a Warner Valve Control device.	MCS_TCX_DEVICE	Connect to an MCS TCX device.
MCS_MEA_SWITCH_DEVICE MCS_MEA_IMPEDANCE_DEVICE MCS_CHANNELTEST_DEVICE MCS_CHANNELTEST_DEVICE MCS_CHANNELTEST_DEVICE MCS_CHANNELTEST_DEVICE MCS_SW2TO64_DEVICE MCS_RETINA_AMS_DONGLE MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOLICE MCS_ROBOOLI	MCS_FCX_DEVICE	Connect to an MCS FCX device.
MCS_MEA_IMPEDANCE_DEVICE MCS_CHANNELTEST_DEVICE MCS_CHANNELTEST_DEVICE MCS_SW2TO64_DEVICE Connect to an MCS ChannelTest device. MCS_RETINA_AMS_DONGLE MCS_RETINA_AMS_DONGLE MCS_RETINA_AMS_DONGLE MCS_RETINA_AMS_DONGLE MCS_RATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE Connect to an MCS Retina AMS Dongle (Radio device) MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOCYTE2_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOLYTE2_DEVICE Connect to an MCS Roboocyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Roboolnject device. MCS_HICLAMP_DEVICE Connect to an MCS Roboolnject device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_MEASURETABLE_DEVICE Connect to an MCS FYI device. MCS_MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_PRS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_D_TESTER_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternBCTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an ExternDTester device. MCS_MEA_CLEAN_DEVICE Connect to an MEA Clean device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a Metal BatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE Connect to a Warner Valve Control device.	MCS_RETINA_LED_DEVICE	Connect to an MCS RetineLed device.
MCS_CHANNELTEST_DEVICE MCS_SW2TO64_DEVICE Connect to an MCS SW2TO64 device. MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS Dongle (Radio device) MCS_PATHIDENT_DEVICE MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE MCS_ROBOO_VEVICE Connect to an MCS Robo Platform device. MCS_ROBOOLEVICE MCS_ROBOOLEVICE MCS_ROBOOLEVICE Connect to an MCS Robooled device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robooled device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robooled device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolly device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolly device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolly device. MCS_ROBOINJECT_DEVICE Connect to an MCS HiClamp device. MCS_HICLAMP_DEVICE Connect to an MCS HiClamp device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_FYI_DEVICE Connect to an MCS Encapsulator device. MCS_HLA_DEVICE Connect to an MCS HLA device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS Peristaltic Pump device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS Peristaltic Pump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_DTESTER_DEVICE Connect to an ExternBCTester device. MCS_MCS_MEA_CLEAN_DEVICE Connect to an MEA Clean device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MBCA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBCA_COAT_DEVICE Connect to a MEA Clean device. MCS_PEC_DOTER_DEVICE Connect to a Mean Clean device. MCS_PEC_DOTER_DEVICE Connect to a Mean Clean device. MCS_PEC_DOTER_DEVICE Connect to a Mean Clean device. MCS_PEC_DOTER_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE Connect to a Warner Valve Control de		Connect to an MCS Mea Switch device.
MCS_SW2TO64_DEVICE MCS_RETINA_AMS_DONGLE Connect to an MCS Retina AMS Dongle (Radio device) MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOO_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOINJECT_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOINJECT_DEVICE Connect to an MCS Roboolyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robologyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robologyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS HiClamp device. MCS_PATCHSERVER_DEVICE Connect to an MCS HiClamp device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_MCS_FYI_DEVICE Connect to an MCS Encapsulator device. MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_FYI_DEVICE Connect to an MCS PYS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS Peristaltic Pump device. MCS_SAFEIS_DEVICE Connect to an MCS Peristaltic Pump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternBCTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an ExternBTester device. MCS_MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MBAEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBAEA_COAT_DEVICE Connect to a MEA Clean device. MCS_PEDOTER_DEVICE Connect to a MaltiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE Connect to a Warner Valve Control device.	MCS_MEA_IMPEDANCE_DEVICE	Connect to an MCS Mea Impedance device.
MCS_RETINA_AMS_DONGLE MCS_PATHIDENT_DEVICE MCS_ROBO_DEVICE MCS_ROBO_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOOCYTE2_DEVICE MCS_ROBOINJECT_DEVICE MCS_ROBOINJECT_DEVICE MCS_ROBOINJECT_DEVICE MCS_ROBOINJECT_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_FYI_DEVICE MCS_FYI_DEVICE MCS_FYI_DEVICE MCS_FYI_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_OKUVISION_STIMULATOR_DEVICE MCS_PREISTALTIC_PUMP_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_CLEAN_DEVICE Connect to an ExternBCTester device. MCS_MEA_CLEAN_DEVICE MCS_MARTIMPLANT_DEVICE MCS_MARTIMPLANT_DEVICE MCS_MARNER_USSING_DEVICE Connect to a PPC device. MCS_PPC_DEVICE Connect to a PPC device. MCS_PPC_DEVICE Connect to a PPC device. MCS_PPC_DEVICE COnnect to a PPC device. MCS_PCCONTROL_DEVICE CONNECT TO A MEA Clean device. MCS_MARNER_USSING_DEVICE CONNECT TO A Warner Valve Control device. MCS_PPC_DEVICE CONNECT TO A MEA Clean device. MCS_MARNER_USSING_DEVICE CONNECT TO A Warner Valve Control device. CONNECT TO A WARNER VALVE_CONTROL_DEVICE CONNECT TO A Warner Valve Control device. CONNECT TO A WARNER VALVE_CONTROL_DEVICE CONNECT TO A WARNER VALVE CONTROL DEVICE CONNECT TO A MARNER VALVE CONTROL DEVICE CONNECT TO A MCS PATITURE. TO A MCS_PROLETAR DEVICE CONNECT TO A MCS PROLETAR CONNEC	MCS_CHANNELTEST_DEVICE	Connect to an MCS ChannelTest device.
MCS_PATHIDENT_DEVICE Connect to an MCS Pathident device. MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOCYTE2_DEVICE Connect to an MCS Roboocyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolnject device. MCS_HICLAMP_DEVICE Connect to an MCS HiClamp device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS PatchServer device. MCS_MEASURETABLE_DEVICE Connect to an MCS Encapsulator device. MCS_FY_DEVICE Connect to an MCS FYI device. MCS_FY_DEVICE Connect to an MCS PYI device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_MCS_MEA_CLEAN_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MARTIMPLANT_DEVICE Connect to a MEA Clean device. MCS_MCS_PEDOTER_DEVICE Connect to a Mea Clean device. MCS_PEDOTER_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PEDOTER_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device.	MCS_SW2TO64_DEVICE	Connect to an MCS SW2TO64 device.
MCS_ROBO_DEVICE Connect to an MCS Robo Platform device. MCS_ROBOOCYTE2_DEVICE Connect to an MCS Roboocyte2 device. MCS_ROBOINJECT_DEVICE Connect to an MCS Robolnject device. MCS_HICLAMP_DEVICE Connect to an MCS HiClamp device. MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_MEASURETABLE_DEVICE Connect to an MCS Encapsulator device. MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_HLA_DEVICE Connect to an MCS PYI device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MBA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBCS_MBA_COAT_DEVICE Connect to a MILIBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PPC device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device.	MCS_RETINA_AMS_DONGLE	Connect to an MCS Retina AMS Dongle (Radio device)
MCS_ROBOCCYTE2_DEVICE	MCS_PATHIDENT_DEVICE	Connect to an MCS Pathident device.
MCS_ROBOINJECT_DEVICE MCS_HICLAMP_DEVICE MCS_PATCHSERVER_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_MEASURETABLE_DEVICE MCS_HILA_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_OKUVISION_STIMULATOR_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_COAT_DEVICE MCS_MEA_COAT_DEVICE MCS_MEA_COAT_DEVICE MCS_MEA_COAT_DEVICE MCS_PPC_DEVICE MCS_PRED_DEVICE MCS_MEA_COAT_DEVICE MCS_PPC_DEVICE MCS_PPC_DEVICE MCS_PROMED To a MEA Clean device. MCS_PPC_DEVICE MCS_PPC_ODEVICE MCS_PPC_DEVICE MCS_PPC_DEVICE MCS_PPC_DEVICE MCS_MEA_CONTROL_DEVICE MC	MCS_ROBO_DEVICE	Connect to an MCS Robo Platform device.
MCS_HICLAMP_DEVICE MCS_PATCHSERVER_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_MEASURETABLE_DEVICE MCS_HLA_DEVICE MCS_FYI_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_PPS_DEVICE MCS_NF_GEN_DEVICE MCS_SAFEIS_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEACOAT_DEVICE MCS_MECS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_MECOS_DEVICE MCS_PECOTER_DEVICE Connect to a MEA Clean device. MCS_MECOS_DEVICE Connect to a MEA Clean device. MCS_MECOS_DEVICE Connect to a MEA Clean device. MCS_MECOS_DEVICE Connect to a Mean Clean device. Connect to a Mean Clean device. MCS_MECOS_DEVICE Connect to a Mean Clean device. MCS_PECOTER_DEVICE Connect to a MultiBatteryCharger device. MCS_PECOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE Connect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device.	MCS_ROBOOCYTE2_DEVICE	Connect to an MCS Roboocyte2 device.
MCS_PATCHSERVER_DEVICE Connect to an MCS PatchServer device. MCS_ENCAPSULATOR_DEVICE Connect to an MCS Encapsulator device. MCS_MEASURETABLE_DEVICE Connect to an MCS Encapsulator device. MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_HLA_DEVICE Connect to an MCS FYI device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS PeristalticPump device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an MEA Clean device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBCOS_DEVICE Connect to a MEA Clean device. MCS_MBCOS_DEVICE Connect to a MILIBAITURED device. MCS_MBCOS_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device.	MCS_ROBOINJECT_DEVICE	Connect to an MCS Robolnject device.
MCS_ENCAPSULATOR_DEVICE MCS_MEASURETABLE_DEVICE MCS_MEASURETABLE_DEVICE MCS_FYI_DEVICE MCS_FYI_DEVICE MCS_HLA_DEVICE MCS_HLA_DEVICE MCS_PPS_DEVICE MCS_NF_GEN_DEVICE MCS_NF_GEN_DEVICE MCS_NF_GEN_DEVICE MCS_SAFEIS_DEVICE MCS_PERISTALTIC_PUMP_DEVICE MCS_EXTERN_BC_TESTER_DEVICE MCS_EXTERN_D_TESTER_DEVICE MCS_SOFTWARE_DONGLE_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_CLEAN_DEVICE MCS_MEA_COAT_DEVICE MCS_MBCOB_DEVICE MCS_MBCOB_DEVICE MCS_MBCOB_DEVICE MCS_MBCOB_DEVICE MCS_MBCOB_DEVICE MCS_MBCOB_DEVICE MCS_PEDOTER_DEVICE Connect to a MEA Clean device. MCS_MBCOB_DEVICE Connect to a SmartImplant device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_HICLAMP_DEVICE	Connect to an MCS HiClamp device.
MCS_MEASURETABLE_DEVICE Connect to an MCS Encapsulator device. MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_HLA_DEVICE Connect to an MCS HLA device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS_DEVICE Connect to an MCS PPS5 device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an Okuvision Stimulator device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafeIS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBCOB_DEVICE Connect to a MultiBatteryCharger device. MCS_MBCOB_DEVICE Connect to a Pedoter device. MCS_PEDOTER_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PATCHSERVER_DEVICE	Connect to an MCS PatchServer device.
MCS_FYI_DEVICE Connect to an MCS FYI device. MCS_HLA_DEVICE Connect to an MCS HLA device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS5_DEVICE Connect to an MCS PPS5 device. MCS_PPS5_DEVICE Connect to an MCS PPS5 device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS PeristalticPump device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBC08_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_ENCAPSULATOR_DEVICE	Connect to an MCS Encapsulator device.
MCS_HLA_DEVICE Connect to an MCS HLA device. MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS5_DEVICE Connect to an MCS PPS5 device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an MCS NF-Gen device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MBCOS_DEVICE Connect to a SmartImplant device. MCS_MBCOS_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a Warner Valve Control device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device.	MCS_MEASURETABLE_DEVICE	Connect to an MCS Encapsulator device.
MCS_PPS_DEVICE Connect to an MCS PPS device. MCS_PPS5_DEVICE Connect to an MCS PPS5 device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an Okuvision Stimulator device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PEDOTER_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_FYI_DEVICE	Connect to an MCS FYI device.
MCS_PPS5_DEVICE Connect to an MCS PPS5 device. MCS_OKUVISION_STIMULATOR_DEVICE Connect to an Okuvision Stimulator device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBCO8_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_HLA_DEVICE	Connect to an MCS HLA device.
MCS_OKUVISION_STIMULATOR_DEVICE Connect to an Okuvision Stimulator device. MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafeIS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PPS_DEVICE	Connect to an MCS PPS device.
MCS_NF_GEN_DEVICE Connect to an MCS NF-Gen device. MCS_SAFEIS_DEVICE Connect to an MCS SafelS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PPS5_DEVICE	Connect to an MCS PPS5 device.
MCS_SAFEIS_DEVICE Connect to an MCS SafeIS device. MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PEDOTER_DEVICE COnnect to a Pedoter device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_OKUVISION_STIMULATOR_DEVICE	Connect to an Okuvision Stimulator device.
MCS_PERISTALTIC_PUMP_DEVICE Connect to an MCS PeristalticPump device. MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.		Connect to an MCS NF-Gen device.
MCS_EXTERN_BC_TESTER_DEVICE Connect to an ExternBCTester device. MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_SAFEIS_DEVICE	Connect to an MCS SafeIS device.
MCS_EXTERN_D_TESTER_DEVICE Connect to an ExternDTester device. MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PERISTALTIC_PUMP_DEVICE	Connect to an MCS PeristalticPump device.
MCS_SOFTWARE_DONGLE_DEVICE Connect to an Software Dongle device. MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.		Connect to an ExternBCTester device.
MCS_MEA_CLEAN_DEVICE Connect to a MEA Clean device. MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_EXTERN_D_TESTER_DEVICE	Connect to an ExternDTester device.
MCS_MEA_COAT_DEVICE Connect to a MEA Clean device. MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_SOFTWARE_DONGLE_DEVICE	Connect to an Software Dongle device.
MCS_SMARTIMPLANT_DEVICE Connect to a SmartImplant device. MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.		
MCS_MBC08_DEVICE Connect to a MultiBatteryCharger device. MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.		
MCS_PEDOTER_DEVICE Connect to a Pedoter device. MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_SMARTIMPLANT_DEVICE	·
MCS_PPC_DEVICE COnnect to a PPC device. WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_MBC08_DEVICE	Connect to a MultiBatteryCharger device.
WARNER_VALVE_CONTROL_DEVICE COnnect to a Warner Valve Control device. WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PEDOTER_DEVICE	Connect to a Pedoter device.
WARNER_USSING_DEVICE COnnect to a Warner Valve Control device.	MCS_PPC_DEVICE	COnnect to a PPC device.
	WARNER_VALVE_CONTROL_DEVICE	COnnect to a Warner Valve Control device.
HEKA_LIH3_DEVICE Connect to a HEKA LIH3 device.	WARNER_USSING_DEVICE	
	HEKA_LIH3_DEVICE	Connect to a HEKA LIH3 device.

ALA_VC3_DEVICE	Connect to an ALA VC3 Valve Commander.
MCS_DEVICE_USB_CYPRESS	Connect to a Cypress USB device.

$\textbf{10.2.1.22} \quad \textbf{DigitalDatastreamEnableEnumNet} \quad \texttt{enum DigitalDatastreamEnableEnumNet} \quad \texttt{[strong]}$

enumerates the streams available as digital datastream

Enumerator

No digital datastream.
16 bits from the standard MUX datastream.
The 16 bits of the standard MUX datastream used by the other virtual device.
The lower 16 bits of the Digital IN port, these ports are on the device by default used as
Digital OUT, thus not available as Digital IN.
The upper 16 bits of the Digital IN port, use when Digital IN datastream is needed.
The lower 16 bits of the Digital OUT port, use when Digital OUT datastream is needed.
The upper 16 bits of the Digital OUT port, these ports are on the device by default used
as Digital IN, thus not available as Digital OUT.
The lower 16 bits of the digital register.
The upper 16 bits of the digital register.
The lower 16 bits of the realtime feedback register.
The upper 16 bits of the realtime feedback register.
The 2 bits of the AUX port.
The 8 bits of the Periodic Pulse Generator (Video-Sync).
The 16 bits of the Digital Out Stimulator.
Headstage 1 digital signals.
Headstage 1 trigger signals.
Headstage 1 lower 16 bits of sideband data.
Headstage 1 upper 16 bits of sideband data.
Headstage 2 digital signals.
Headstage 2 trigger signals.
Headstage 2 lower 16 bits of sideband data.
Headstage 2 upper 16 bits of sideband data.

10.2.1.23 DigitalStimulatorTriggerEventEnumNet enum DigitalStimulatorTriggerEventEnumNet [strong]

Enumerates start/stop event for DigOut/DigStim trigger. /summary>

Start	
Stop	

10.2.1.24 DigitalStimulatorTriggerSlopeEnumNet enum DigitalStimulatorTriggerSlopeEnumNet [strong]

Enumerates start/stop conditions for DigOut/DigStim trigger. /summary>

Enumerator

Falling	
Rising	

10.2.1.25 DigitalTargetEnumNet enum DigitalTargetEnumNet [strong]

Enumerates the Digital Targets for Digital Sources

Enumerator

10.2.1.26 ElectrodeDacMuxEnumNet enum ElectrodeDacMuxEnumNet [strong]

Enumerates the setting of the Stimulation DAC Multiplexer.

Enumerator

	Ground	Connect the electrode to Ground while stimulation is active.
Stg1 Connect the electrode tp STG 1 while stimulation is ac		Connect the electrode tp STG 1 while stimulation is active.
Stg2 Connect the electrode tp STG 2 while stimulation is active		Connect the electrode tp STG 2 while stimulation is active.
	Stg3	Connect the electrode tp STG 3 while stimulation is active.

10.2.1.27 ElectrodeModeEnumNet enum ElectrodeModeEnumNet [strong]

Enumerates the mode of each electrode, can be automatic or manual. In automatic mode, the blanking of the electrode is controlled by the sideband signal, in manual mode, the stimulation configuration is independant of the sideband signal.

emAutomatic	
emManual	

10.2.1.28 enCMosMeaChipType enum enCMosMeaChipType [strong]

Enumerator

unknown	
nMos16LV	
nMos32LV	
nMos36LN	
nMos64LN	

10.2.1.29 EnSTG200x_STATUS enum EnSTG200x_STATUS [strong]

Enumerator

OK	
NOT_CONNECTED	
DEVICE_NOT_FOUND	

10.2.1.30 FilterAttributeEnumNet enum FilterAttributeEnumNet [strong]

Enumerator

PreCommaB	
PostCommaB	
CommaPositionB	
PreCommaA	
PostCommaA	
CommaPositionA	

10.2.1.31 FilterBandEnumNet enum FilterBandEnumNet [strong]

Enumerator

Unknown

Lowpass	
Highpass	

10.2.1.32 FilterCalculationDirectionEnumNet enum FilterCalculationDirectionEnumNet [strong]

Enumerator

DoubleToInt	
IntToDouble	

10.2.1.33 FilterFamilyEnumNet enum FilterFamilyEnumNet [strong]

Enumerator

Unknown	
Bessel	
Butterworth	
RC	

10.2.1.34 FilterTypeEnumNet enum FilterTypeEnumNet [strong]

Enumerator

Hardware	
Software	

10.2.1.35 FpgaldEnumNet enum FpgaldEnumNet [strong]

DeviceNotConnected	
Mea2100Interfaceboard	
Mea2100Headstage	
Mea2100STG	
MultiwellHeadstage	
MultiwellInterfaceboard	

TbsiDacqInterfaceboard	
TbsiDacqHeadstage	
CmosMeaInterfaceboard	
CmosMeaHeadstage	
Mea2100MultiwellIFB2	
Me2100Interfaceboard	
Me2100InvivoSignalCollectorUnit	
Me2100InvitroSignalCollectorUnit	
Me2100_32XilinxHeadstage	
Me2100_32PICiCE40Headstage	
Mea2100_256Interfaceboard	
Mea2100_256Headstage	
W2100Interfaceboard	
W2100WirelessReceiver	
W2100WirelessReceiverAnalog	
Mea2100Mini60PICiCE40Headstage	
Mea2100BetaScreenHeadstage	
Me2100UPA32Headstage	
MultiwellMiniHeadstage	
Mea2100Mini120Headstage	
Mea2100Mini60ECP5Headstage	
eCubeHeadstage	
Me2100Graphene16_32Headstage	
GrapheneASICHeadstage	
WholeCellPatchHeadstage	
InterfaceBoard2	
W2100IFB2	
CmosmealFB2	
M4M_Dock	
M4M_ReUse	
Mea2100LiteHeadstage	
LIH30Interfaceboard	
LIH30ADCCtrl	
UssingRail	
UssingChamber	
IFB2GoldenInterfaceboard	
IFB30GoldenInterfaceboard	
DeviceHasNoHeadstage	

10.2.1.36 HeadstageldEnumNet enum HeadstageIdEnumNet [strong]

DeviceNotConnected	
Mea2100	

Multiwell	
TbsiDacq	
CmosMea	
InvivoSignalCollectorUnit	
InvitroSignalCollectorUnit	
Mea2100_256	
W2100WirelessReceiver	
W2100WirelessReceiverAnalog	
Mea2100_Lite	
LIH30ADCCtrl	
DeviceHasNoHeadstage	

10.2.1.37 loVoltageEnumNet enum IoVoltageEnumNet [strong]

enumerates the I/O Voltages available on the IFB2

Enumerator

Voltage_3V3	
Voltage_5V0	

10.2.1.38 LIH30_ADC_Channel_EnumNet enum LIH30_ADC_Channel_EnumNet [strong]

User_ADC_0	
User_ADC_1	
User_ADC_2	
User_ADC_3	
User_ADC_4	
Test_ADC_EPC10	
ModulA_ADC0	
ModulA_ADC1	
ModulA_ADC2	
ModulA_ADC3	
ModulB_ADC0	
ModulB_ADC1	
ModulB_ADC2	
ModulB_ADC3	
ModulC_ADC0	
ModulC_ADC1	
ModulC_ADC2	
ModulC_ADC3	

ModulD_ADC0	
ModulD_ADC1	
ModulD_ADC2	
ModulD_ADC3	

10.2.1.39 LIH30_DAC_Channel_EnumNet enum LIH30_DAC_Channel_EnumNet [strong]

Enumerator

User_DAC_0	
User_DAC_1	
User_DAC_2	
Test_DAC_EPC10	
ModulA_DAC0	
ModulA_DAC1	
ModulB_DAC0	
ModulB_DAC1	
ModuIC_DAC0	
ModulC_DAC1	
ModuID_DAC0	
ModuID_DAC1	

10.2.1.40 LIH30_EPC10_Bus_EnumNet enum LIH30_EPC10_Bus_EnumNet [strong]

Enumerator

Α	
В	

10.2.1.41 MbcChannelStateEnumNet enum MbcChannelStateEnumNet [strong]

Enumerator

csIdleNoBattery csIdleChargeFinished csCapacityTestPrecharge csCapacityTestDischarge csRefreshBattery csCharge csDischarge csError		
csCapacityTestPrecharge csCapacityTestDischarge csRefreshBattery csCharge csDischarge	csldleNoBattery	
csCapacityTestDischarge csRefreshBattery csCharge csDischarge	csIdleChargeFinished	
csRefreshBattery csCharge csDischarge	csCapacityTestPrecharge	
csCharge csDischarge	csCapacityTestDischarge	
csDischarge	csRefreshBattery	
	csCharge	
csError	csDischarge	
	csError	

Generated by Doxygen

10.2.1.42 MbcChargingModeEnumNet enum MbcChargingModeEnumNet [strong]

Enumerator

StorageCharge	
FullCharge	

10.2.1.43 MbcRatedCapacityEnumNet enum MbcRatedCapacityEnumNet [strong]

Enumerator

rc30mAh	
rc100mAh	
rc200mAh	
rc300mAh	
rcGreater300mAh	

10.2.1.44 McsBusTypeEnumNet enum McsBusTypeEnumNet [strong]

Enumerates the bus to use, either USB, PCI or any

Enumerator

MCS_ANY_BUS	
MCS_UNDEFINED_BUS	
MCS_USB_BUS	
MCS_PCI_BUS	

10.2.1.45 McsUsbSpeedEnumNet enum McsUsbSpeedEnumNet [strong]

Enumerates the current connection speed of the device

LowSpeed	
FullSpeed	
HighSpeed	
SuperSpeed	
UnknownSpeed	

10.2.1.46 MEA2100_256DacqGroupChannelEnumNet enum MEA2100_256DacqGroupChannelEnumNet [strong]

Enumerates the MEA2100-256 Channel Groups of Datastream

Enumerator

HS1ElectrodeGroup	
HS2ElectrodeGroup	
InterfaceADCGroup	
STG1DACSignalGroup	
STG2DACSignalGroup	
DSPAnalogGroup	
DSPDigitalGroup	
IFDigChannelsGroup	
STG1TriggerStatusGroup	
STG1SidebandsGroup	
STG2TriggerStatusGroup	
STG2SidebandsGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.47 MEA2100_256DigitalSourceEnumNet enum MEA2100_256DigitalSourceEnumNet [strong]

Enumerates the digital source of the MEA2100-256 device.

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
DeviceRunStatus	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Trigger3Status	
HS1Trigger4Status	
HS1Trigger5Status	
HS1Trigger6Status	
HS1Trigger7Status	

Eliumerator	
HS1Trigger8Status	
HS1Trigger9Status	_
HS1Trigger10Status	_
HS1Trigger11Status	_
HS1Trigger12Status	_
HS1Trigger13Status	
HS1Trigger14Status	
HS1Trigger15Status	
HS1Trigger16Status	
HS1Trigger17Status	
HS1Trigger18Status	
HS1Sideband1	_
HS1Sideband2	_
HS1Sideband3	_
HS1Sideband4	_
HS1Sideband5	
HS1Sideband6	_
HS1Sideband7	_
HS1Sideband8	_
HS1Sideband9	_
HS1Sideband10	
HS1Sideband11	_
HS1Sideband12	_
HS1Sideband13	_
HS1Sideband14	
HS1Sideband15	_
HS1Sideband16	
HS1Sideband17	
HS1Sideband18	
HS2Trigger1Status	
HS2Trigger2Status	_
HS2Trigger3Status	_
HS2Trigger4Status	_
HS2Trigger5Status	_
HS2Trigger6Status	_
HS2Trigger7Status	
HS2Trigger8Status	
HS2Trigger9Status	
HS2Trigger10Status	
HS2Trigger11Status	_
HS2Trigger12Status	
HS2Trigger13Status	
HS2Trigger14Status	_
HS2Trigger15Status	_
HS2Trigger16Status	_
HS2Trigger17Status	_
	_
HS2Trigger18Status	
HS2Sideband1	

HS2Sideband2	
HS2Sideband3	
HS2Sideband4	
HS2Sideband5	
HS2Sideband6	
HS2Sideband7	
HS2Sideband8	
HS2Sideband9	
HS2Sideband10	
HS2Sideband11	
HS2Sideband12	
HS2Sideband13	
HS2Sideband14	
HS2Sideband15	
HS2Sideband16	
HS2Sideband17	
HS2Sideband18	
LastPosition	

10.2.1.48 MEA2100DigitalSourceEnumNet enum MEA2100DigitalSourceEnumNet [strong]

Enumerates the digital source of the MEA2100 device.

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Trigger3Status	
HS1Trigger4Status	
HS1Trigger5Status	
HS1Trigger6Status	
HS1Sideband1	
HS1Sideband2	
HS1Sideband3	
HS1Sideband4	
HS1Sideband5	
HS1Sideband6	
HS2Trigger1Status	
HS2Trigger2Status	
HS2Trigger3Status	

HS2Trigger4Status	
HS2Trigger5Status	
HS2Trigger6Status	
HS2Sideband1	
HS2Sideband2	
HS2Sideband3	
HS2Sideband4	
HS2Sideband5	
HS2Sideband6	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
LastPosition	

Enumerates the MEA layout of the MEA2100 device.

Enumerator

mlUnknown	
mIMEA60	

10.2.1.50 MultiwellDigitalSourceEnumNet enum MultiwellDigitalSourceEnumNet [strong]

Enumerates the digital source of the Multiwell device.

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
HS1Trigger1Status	
HS1Trigger2Status	
HS1Sideband1	
HS1Sideband2	
HS2Trigger1Status	
HS2Trigger2Status	

HS2Sideband1	
HS2Sideband2	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	
LastPosition	

10.2.1.51 MultiwellPlateTypeEnumNet enum MultiwellPlateTypeEnumNet [strong]

Enumerator

Plate_Dummy	
Plate_24W700_100FMA	
Plate_24W030MGA	
Plate_72W500_100PMA	
Plate_72W500_100FMA	
Plate_24W700_100FMB	
Plate_96W700_100FMA	
Plate_96W300_80_1152FMA	
Plate_96W400_80_1152FMB	
Plate_24W300_30_1152GBA	
Plate_24W700_100FMC	
Plate_96W700_100FMB	
Plate_96W700_100GBC	
Plate_96W700_100GBD	
Plate_24W700_100PBA	
Plate_Dummy_126	
Plate_24W300_30GMA	
Plate_96W700_100GMA	
Plate_24W300_30GBA	
Plate_96W700_100GBA	
Plate_24W300_30GBB	
Plate_96W700_100GBB	
No_Plate	

10.2.1.52 PatchServAdcModeEnumNet enum PatchServAdcModeEnumNet [strong]

Normal	
CatchAmp	

10.2.1.53 PlateClampEnumNet enum PlateClampEnumNet [strong]

Enumerator

Close	
Open	
Stop	

10.2.1.54 PlateClampLockEnumNet enum PlateClampLockEnumNet [strong]

Enumerator

Lock	
Unlock	

10.2.1.55 PortDirectionEnumNet enum PortDirectionEnumNet [strong]

enumerates a port direction

Enumerator

Output	
Input	

10.2.1.56 PP_Pump_Mode_Type_EnumNet enum PP_Pump_Mode_Type_EnumNet [strong]

Enumerator

Manual	
Digital	
Analog	

10.2.1.57 ProductIdEnumNet enum ProductIdEnumNet [strong]

Enumerates the group of MCS devices to connect to.

Any	
None	
LegacyMeaUsb	
ALA_VC3	
Cypress_FX1	
Cypress_FX2	
Cypress_FX3	
MC Card	
Campden_Ci4600EphysVideoDataIntegrator	
HekaLIH30	
HekaEPC10Single	
HekaEPC10Double	
HekaEPC10Triple	-
HekaEPC10Quadro	-
HekaLIH406	-
HekaLIH816	
HekalTEV100	
HekaPG610	
HekaPG611	-
HekaPG612	
HekaPG618	
HekaPG690	
HekaEPCLite	
STG	
Octopot	
Tersens	
Dotriapot	
HLA	
STG400x	
STG4002	
STG4004	
STG4008	
STG400x_opto	
STG4002_opto	
STG4004_opto	
STG4008_opto	
STG5	
STG3008_FA	
MultiwellOptoStim	
Generic	
PGA	
PCX	
TCX	
FCX	
FCB	
TC01	
TC02	
Retina_LED	
AMS_Dongle	

Okuvision Stimulator	
ExternBCTester	
Triggerbox_IMS	
Triggerbox_AMS	
Triggerbox_AMS3	
ExternDTester	
FunkDongleS	
ExternSTester	
DongleS	
Triggerbox_R5	
MEA Switch	
MEA_Impedance	
ChannelTest	
Sw2to64	
PeristalticPump	
MEA Switch 2 1	
MEA_Switch_2_1 MEA Switch 4 2	
PPS4plus1	
PPS5	
PPS2	
PPS5_DIG	
MEA_Clean	
MEA_Coat	
Multiwell_ICC	
MBC08	
PPC	
MEA1060	
MEA_Sanofi	
ME256	
ME128	
ME64	
ME32	
ME16	
MEA2100_Mini_Usb_develop	
MEA256	
MEA2100	
MEA2100_32	
MEA2100_Lite	
Multiwell	
MEA2100_256	
ME2100	
MEA2100BetaScreen	
MEA2100_Mini	
TBSI_Dacq	
Multiwell MEA Mini	
Whole Cell Patch	
eCube	
Graphene_ASIC	
GE2100	
GL2100	

Multiboot	
WPA8	
WPA4	
WPA16	
WPA32	
W2100	
NeuroChip	
UsbTest	
SoftwareDongle	
PathIdent	
NF_Gen	
SafelS	
Encapsulator	
NeurochipConfig	
MeasureTable	
Robooycte2	
Robolnject	
HiClamp	
PatchServer	
Dilutor	
HiClamp4Uart	
IM16S16KRA	
IM64KRB	
IS32KRA	
IM64KRC	
IM16S8KRA	
IM16KRC	
SmartImplant	
PositionImp	
PositionBase	
PositionIICentralUnit	
PositionIIBase	
GrapheneProjectTestDevice	
Pos900	
Neptun	
Warner_Valve_Control	
Warner_TEER_Machine	
Warner_Ussing	

10.2.1.58 PulseGenerator_Mode_EnumNet enum PulseGenerator_Mode_EnumNet [strong]

Off	
AlwaysOn	
Gated_Low_Active	
Gated_High_Active	

10.2.1.59 ReferenceElectrodeModeEnumNet enum ReferenceElectrodeModeEnumNet [strong]

enumerates the electrode subtraction modes

Enumerator

SubtractionOff	
SubtractFromAllOther	
SubtractFromReferenceElectrodeOnly	
SubtractFromAll	

10.2.1.60 ReferenceElectrodeSwitchPositionEnumNet enum ReferenceElectrodeSwitchPositionEnumNet [strong]

enumerates the possible positions of the reference electrode switch of the ME2100 device

Enumerator

off	
Ref8	
Ref16	
Ref24	
Ref32	

10.2.1.61 RetriggerActionEnumNet enum RetriggerActionEnumNet [strong]

Enumerates possible retrigger actions for STG200x devices.

Enumerator

Stop	
Restart	
Ignore	
Gate	
Single	

10.2.1.62 RoboCurrentModeEnumNet enum RoboCurrentModeEnumNet [strong]

Off	
Break	
Standby	
Reference	
Movement	

10.2.1.63 SampleDstSizeNet enum SampleDstSizeNet [strong]

Enumerates the destination data format for ChannelBlock functions.

Enumerator

SampleDstSize16	
SampleDstSize32	

10.2.1.64 SampleSizeNet enum SampleSizeNet [strong]

Enumerates the data format for ChannelBlock functions.

Enumerator

SampleSize16Unsigned	
SampleSize16Signed	
SampleSize24Unsigned	
SampleSize24Signed	
SampleSize32Unsigned	
SampleSize32Signed	
SampleSize64Unsigned	
SampleSize64Signed	

10.2.1.65 SCU_HeadstageldEnumNet enum SCU_HeadstageIdEnumNet [strong]

DeviceNotConnected	
Me2100_32Xilinx	
Me2100_32PICiCE40	
Mea2100Mini60PICiCE40	
Mea2100BetaScreen	
Me2100UPA32	

MultiwellMini	
Mea2100Mini120	
Mea2100Mini60ECP5	
eCube	
Me2100Graphene16_32	
GrapheneASIC	
WholeCellPatch	
DeviceHasNoHeadstage	

10.2.1.66 SCUDacqGroupChannelEnumNet enum SCUDacqGroupChannelEnumNet [strong]

Enumerates the SCU Channel Groups of Datastream

Enumerator

10.2.1.67 SCUDigitalSourceEnumNet enum SCUDigitalSourceEnumNet [strong]

Enumerates the digital source of the SCU device.

Enumerator	
DigitalInOfOutPort	
DigitalIn	_
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
DeviceRunStatus	_
SCU1HS1Sideband1TriggerStatus	_
SCU1HS1Sideband2TriggerStatus	_
SCU1HS2Sideband1TriggerStatus	_
SCU1HS2Sideband2TriggerStatus	
SCU1HS3Sideband1TriggerStatus	_
SCU1HS3Sideband2TriggerStatus	
SCU1HS4Sideband1TriggerStatus	
SCU1HS4Sideband2TriggerStatus	
SCU1LEDStim1TriggerStatus	
SCU1LEDStim2TriggerStatus	
SCU1LEDStim3TriggerStatus	
SCU1LEDStim4TriggerStatus	
SCU1HS1Sideband1	
SCU1HS1Sideband2	
SCU1HS2Sideband1	
SCU1HS2Sideband2	
SCU1HS3Sideband1	
SCU1HS3Sideband2	
SCU1HS4Sideband1	
SCU1HS4Sideband2	
SCU1LEDStim1	
SCU1LEDStim2	
SCU1LEDStim3	
SCU1LEDStim4	
SCU2HS1Sideband1TriggerStatus	
SCU2HS1Sideband2TriggerStatus	
SCU2HS2Sideband1TriggerStatus	_
SCU2HS2Sideband2TriggerStatus	
SCU2HS3Sideband1TriggerStatus	
SCU2HS3Sideband2TriggerStatus	_
SCU2HS4Sideband1TriggerStatus	_
SCU2HS4Sideband2TriggerStatus	
SCU2LEDStim1TriggerStatus	_
SCU2LEDStim2TriggerStatus	_
SCU2LEDStim3TriggerStatus	
SCU2LEDStim4TriggerStatus	
SCU2HS1Sideband1	

SCU2HS1Sideband2	
SCU2HS2Sideband1	
SCU2HS2Sideband2	
SCU2HS3Sideband1	
SCU2HS3Sideband2	
SCU2HS4Sideband1	
SCU2HS4Sideband2	
SCU2LEDStim1	
SCU2LEDStim2	
SCU2LEDStim3	
SCU2LEDStim4	
LastPosition	
	SCU2HS2Sideband1 SCU2HS2Sideband2 SCU2HS3Sideband2 SCU2HS3Sideband2 SCU2HS4Sideband1 SCU2HS4Sideband2 SCU2HS4Sideband2 SCU2LEDStim1 SCU2LEDStim2 SCU2LEDStim3 SCU2LEDStim4

10.2.1.68 Stg200xDigoutModeEnumNet enum Stg200xDigoutModeEnumNet [strong]

Enumerates the DigoutMode on STG400x devices.

Enumerator

Monitor	Monitor digital input pins. Digital out is a mirror of the digital input pins.
Manual	Manually set the value on the digital out pins with SetDigoutValue.
SYNCOUT1	show bit 7 to 15 of syncout channel 1 on the digital outputs
SYNCOUT2	show bit 7 to 15 of syncout channel 2 on the digital outputs
SYNCOUT3	show bit 7 to 15 of syncout channel 3 on the digital outputs
SYNCOUT4	show bit 7 to 15 of syncout channel 4 on the digital outputs
SYNCOUT5	show bit 7 to 15 of syncout channel 5 on the digital outputs
SYNCOUT6	show bit 7 to 15 of syncout channel 6 on the digital outputs
SYNCOUT7	show bit 7 to 15 of syncout channel 7 on the digital outputs
SYNCOUT8	show bit 7 to 15 of syncout channel 8 on the digital outputs

$\textbf{10.2.1.69} \quad \textbf{Stg200xMultiFileSubmodeEnumNet} \quad \texttt{enum Stg200xMultiFileSubmodeEnumNet} \quad \texttt{[strong]}$

Enumerates EnableMultifileMode submodes for STG devices.

Enumerator

MultiFileMode	No Flags.
ExtendedMultiFileMode	Assign all channels to the trigger which number is the given segment number.

$\textbf{10.2.1.70} \quad \textbf{Stg200xSegmentFlagsEnumNet} \quad \texttt{enum Stg200xSegmentFlagsEnumNet} \quad \texttt{[strong]}$

Enumerates Segmentflag options for STG400x devices.

Enumerator

None	No Flags.
UpdateTrigger	Assign all channels to the trigger which number is the given segment number.
DownloadOnly	Only switch the segment for the next download, keep current segment running.
TriggerOnly	Only switch the segment for the next sweep, keep current downlaod segment.
SyncStart	Delay the start the new segment with SegmentStart() until the next sweep has finished.

10.2.1.71 Stg200xTriggerStatusEnumNet enum Stg200xTriggerStatusEnumNet [strong]

Enumerates the STG download mode trigger status

The STG maintains the status for each of the STG200x_NUM_TRIGGER triggers

Enumerator

ldle	
Running	
Finished	
Armed	

$\textbf{10.2.1.72} \quad \textbf{Stg3008FilterAmpAmplificationEnumNet} \quad \texttt{enum Stg3008FilterAmpAmplificationEnumNet} \quad \texttt{[strong]}$

Enumerates the STG3008FA filter amplifier gains

Enumerator

Gain1000	
Gain500	
Gain200	
Gain100	

10.2.1.73 STG_DestinationEnumNet enum STG_DestinationEnumNet [strong]

Enumerates the destination for STG downloads.

channeldata_voltag	е
channeldata_curre	nt
syncoutda	а

channeldata_positive_voltage	
channeldata_positive_current	
rawdata	
channeldata_current_own_sync	
channeldata_positive_current_own_sync	
channeldata_current_own_boost_gnd_sync	
channeldata_positive_current_own_boost_gnd_sync	
channeldata_current_always_boost	
channeldata_current_always_boost_own_sync	

10.2.1.74 StimulationLayoutConfigurationEnumNet enum StimulationLayoutConfigurationEnumNet [strong]

enumerates the layout configuration for the MEA2100-256 device

Enumerator

SingleWell	
SixWell	
NineWell	

10.2.1.75 TBSI_DACQDigitalSourceEnumNet enum TBSI_DACQDigitalSourceEnumNet [strong]

Enumerates the digital source of the TBSI-DACQ device.

DigitalInOfOutPort	
DigitalIn	
DigitalPulse	
Feedback	
AuxIn	
Zero	
One	
DeviceRunStatus	
PulseGenerator	
DigitalOutStimulator	
DigitalData	
HS1DigitalData1	
HS2DigitalData1	
LastPosition	

Enumerates the type of TCX devices.

Enumerator

Unknown	
Regular	
BMI	
Nanion	
Warner	

10.2.1.77 TcxSensorTypeEnumNet enum TcxSensorTypeEnumNet [strong]

Enumerates the sensor types for TCX devices

Enumerator

Reserved5	
Reserved4	
Reserved3	
Reserved2	
Reserved1	
NTC10K	
PT1000	
PT100	

10.2.1.78 TeerClampModeEnumNet enum TeerClampModeEnumNet [strong]

Enumerator

ClampModeVoltage	
ClampModeCurrent	
ClampModeOpen	
ClampModeInternalCalibration	

10.2.1.79 TeerWaveformEnumNet enum TeerWaveformEnumNet [strong]

Rectangle	
Sine	

10.2.1.80 TriggerSourceEnumNet enum TriggerSourceEnumNet [strong]

Enumerates the trigger source of the MEA2100 device.

tsNone tsDigitalIn1 tsDigitalIn2 tsDigitalIn3 tsDigitalIn4 tsDigitalIn5 tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn15 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn19 tsDigitalIn19 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2 tsFeedback3	Litatificiatoi	
tsDigitalIn2 tsDigitalIn3 tsDigitalIn4 tsDigitalIn5 tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn15 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsNone	
tsDigitalIn3 tsDigitalIn4 tsDigitalIn5 tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn1	
tsDigitalIn4 tsDigitalIn5 tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn15 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn2	
tsDigitalIn5 tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn15 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn3	
tsDigitalIn6 tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn4	
tsDigitalIn7 tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn22 tsDigitalIn23 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn5	
tsDigitalIn8 tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsPeedback1 tsFeedback2	tsDigitalIn6	
tsDigitalIn9 tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsPeedback1 tsFeedback2	tsDigitalIn7	
tsDigitalIn10 tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn21 tsDigitalIn22 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn8	
tsDigitalIn11 tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsDigitalIn32 tsDigitalIn32	tsDigitalIn9	
tsDigitalIn12 tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn23 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn10	
tsDigitalIn13 tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsDigitalIn32	tsDigitalIn11	
tsDigitalIn14 tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn25 tsDigitalIn27 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsPeedback1 tsFeedback2	tsDigitalIn12	
tsDigitalIn15 tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsPeedback1 tsFeedback2	tsDigitalIn13	
tsDigitalIn16 tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsPeedback1 tsFeedback2	tsDigitalIn14	
tsDigitalIn17 tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn15	
tsDigitalIn18 tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn16	
tsDigitalIn19 tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn17	
tsDigitalIn20 tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn18	
tsDigitalIn21 tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn19	
tsDigitalIn22 tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn20	
tsDigitalIn23 tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn21	
tsDigitalIn24 tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn22	
tsDigitalIn25 tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn23	
tsDigitalIn26 tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn24	
tsDigitalIn27 tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn25	
tsDigitalIn28 tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn26	
tsDigitalIn29 tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn27	
tsDigitalIn30 tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn28	
tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn29	
tsDigitalIn31 tsDigitalIn32 tsFeedback1 tsFeedback2	tsDigitalIn30	
tsFeedback1 tsFeedback2		
tsFeedback2	tsDigitalIn32	
	tsFeedback1	
tsFeedback3	tsFeedback2	
	tsFeedback3	

tsFeedbac	k4
tsFeedbac	
tsFeedbac	k6
tsFeedbac	k7
tsFeedbac	
tsFeedbac	k9
tsFeedback	10
tsFeedback	11
tsFeedback	12
tsFeedback	
tsFeedback	
tsFeedback	15
tsFeedback	
tsFeedback:	
tsFeedback:	
tsFeedback	
tsFeedback:	
tsFeedback	
tsFeedback	
tsFeedback:	
tsFeedback:	
tsAuxli	
tsAuxli	
tsDigitalPuse	
tsDigitalPuse	
tsDigitalPuse	e2
tsDigitalPuse	e3
tsDigitalPuse	e 4
tsDigitalPuse	e5
tsDigitalPuse	
tsDigitalPuse	
tsDigitalPuse	
tsDigitalPuse	
tsDigitalPuse	14
tsDigitalPuse ²	15
tsDigitalPuse ²	16
tsDigitalPuse	17
<u> </u>	

tsDigitalPuse18	
tsDigitalPuse19	
tsDigitalPuse20	
tsDigitalPuse21	
tsDigitalPuse22	
tsDigitalPuse23	
tsDigitalPuse24	
tsDigitalPuse25	
tsDigitalPuse26	
tsDigitalPuse27	
tsDigitalPuse28	
tsDigitalPuse29	
tsDigitalPuse30	
tsDigitalPuse31	
tsTriggered	
tsSidebandBit8	
tsDACQCy1Dev1Runs	
tsDACQCy1Dev2Runs	
tsDACQCy2Dev1Runs	
tsDACQCy2Dev2Runs	

10.2.1.81 UsbVendorldEnumNet enum UsbVendorldEnumNet [strong]

Enumerator

Unknown	
None	
Renesas	
ASMedia	
Intel	

$\textbf{10.2.1.82} \quad \textbf{UssingClampModeEnumNet} \quad \texttt{enum UssingClampModeEnumNet} \quad \texttt{[strong]}$

VoltageClamp	
CurrentClamp	
OpenClamp	
Standby	
ElectrodeOffset	

10.2.1.83 UssingUnitEnumNet enum UssingUnitEnumNet [strong]

Enumerator

Volt	
Ampere	
State	

10.2.1.84 VendorldEnumNet enum VendorldEnumNet [strong]

Enumerates the group of MCS devices to connect to.

Enumerator

Any	
None	
MCS	
PCI	
Cypress	
ALA_VC3	

10.2.1.85 W2100_Accel_Gyro_Select_EnumNet enum W2100_Accel_Gyro_Select_EnumNet [strong]

enumerates the accelerometer configuration on the W2100 device

Enumerator

Off	
GyroOnly	
AccelOnly	
Both	

10.2.1.86 W2100DacqGroupChannelEnumNet enum W2100DacqGroupChannelEnumNet [strong]

Enumerates the W2100 Channel Groups of Datastream

Enumerator

InterfaceADCGroup

DSPDataGroup	
WirelessHeadStageAnalogRE1HS1	
WirelessHeadStageStatusRE1HS1	
WirelessHeadStageAnalogRE1HS2	
WirelessHeadStageStatusRE1HS2	
WirelessHeadStageAnalogRE1HS3	
WirelessHeadStageStatusRE1HS3	
WirelessHeadStageAnalogRE1HS4	
WirelessHeadStageStatusRE1HS4	
WirelessHeadStageAnalogRE2HS1	
WirelessHeadStageStatusRE2HS1	
WirelessHeadStageAnalogRE2HS2	
WirelessHeadStageStatusRE2HS2	
WirelessHeadStageAnalogRE2HS3	
WirelessHeadStageStatusRE2HS3	
WirelessHeadStageAnalogRE2HS4	
WirelessHeadStageStatusRE2HS4	
WirelessHeadStageGyroDataRE1HS1	
WirelessHeadStageAccDataRE1HS1	
WirelessHeadStageGyroDataRE1HS2	
Wireless HeadStageAccDataRE1HS2	
Wireless HeadStageGyroDataRE1HS3	
WirelessHeadStageAccDataRE1HS3	
Wireless HeadStageGyroDataRE1HS4	
WirelessHeadStageAccDataRE1HS4	
WirelessHeadStageGyroDataRE2HS1	
- 1	
WirelessHeadStageAccDataRE2HS1	
WirelessHeadStageGyroDataRE2HS2	
WirelessHeadStageAccDataRE2HS2	
WirelessHeadStageGyroDataRE2HS3	
WirelessHeadStageAccDataRE2HS3	
WirelessHeadStageGyroDataRE2HS4	
WirelessHeadStageAccDataRE2HS4	
WirelessHeadStageOptoStimCurrentRE1HS1	
WirelessHeadStageReservedARE1HS1	
WirelessHeadStageOptoStimCurrentRE1HS2	
WirelessHeadStageReservedARE1HS2	
WirelessHeadStageOptoStimCurrentRE1HS3	
WirelessHeadStageReservedARE1HS3	
WirelessHeadStageOptoStimCurrentRE1HS4	
WirelessHeadStageReservedARE1HS4	
WirelessHeadStageOptoStimCurrentRE2HS1	
WirelessHeadStageReservedARE2HS1	
WirelessHeadStageOptoStimCurrentRE2HS2	
WirelessHeadStageReservedARE2HS2	
WirelessHeadStageOptoStimCurrentRE2HS3	

WirelessHeadStageReservedARE2HS3	
WirelessHeadStageOptoStimCurrentRE2HS4	
WirelessHeadStageReservedARE2HS4	
WirelessHeadStageReservedBRE1HS1	
WirelessHeadStageReservedCRE1HS1	
WirelessHeadStageReservedBRE1HS2	
WirelessHeadStageReservedCRE1HS2	
WirelessHeadStageReservedBRE1HS3	
WirelessHeadStageReservedCRE1HS3	
WirelessHeadStageReservedBRE1HS4	
WirelessHeadStageReservedCRE1HS4	
WirelessHeadStageReservedBRE2HS1	
WirelessHeadStageReservedCRE2HS1	
WirelessHeadStageReservedBRE2HS2	
WirelessHeadStageReservedCRE2HS2	
WirelessHeadStageReservedBRE2HS3	
WirelessHeadStageReservedCRE2HS3	
WirelessHeadStageReservedBRE2HS4	
WirelessHeadStageReservedCRE2HS4	
IFDigChannelsGroup	
AudioTestChannelGroup	
PacketFrameContextGroup	

10.2.1.87 W2100DigitalSourceEnumNet enum W2100DigitalSourceEnumNet [strong]

Enumerates the digital source of the W2100 device.

10.2.1.88 WvcDisplayModeEnumNet enum WvcDisplayModeEnumNet [strong]

enumerates Wvc display mode

Enumerator

Work	
PC	
Settings	
TouchTest	

10.2.1.89 WvcValveModeEnumNet enum WvcValveModeEnumNet [strong]

enumerates Wvc valve mode

Enumerator

Manual	
Digital	
Analog	
Table	

10.2.2 Function Documentation

```
10.2.2.2 OnDeviceArrivalRemoval() public delegate void Mcs::Usb::OnDeviceArrivalRemoval (
CMcsUsbListEntryNet^ entry)
```

Delegate to show a device arrival or removal.

```
10.2.2.4 OnMcsUsbDeviceState() public delegate void OnMcsUsbDeviceState (
             usbSetupPacket_t^ request )
10.2.2.5 OnMcsUsbDeviceStateCallback() private delegate void OnMcsUsbDeviceStateCallback (
             IntPtr pThis,
             uint32_t size,
             IntPtr buffer )
10.2.2.6 OnMwPollStatus() public delegate void Mcs::Usb::OnMwPollStatus (
             unsigned int CurrentTemp,
             unsigned int PlateState,
             unsigned int SwitchState )
10.2.2.7 OnStg200xDataHandler() public delegate void Mcs::Usb::OnStg200xDataHandler (
             uint32_t trigger )
10.2.2.8 OnStg200xErrorHandler() public delegate void Mcs::Usb::OnStg200xErrorHandler ( )
10.2.2.9 OnStgPollStatus() public delegate void Mcs::Usb::OnStgPollStatus (
             unsigned int status,
             StgStatusNet^ stgStatusNet,
             array < int >^{\land} index_list)
10.2.2.10 OnUpdateFirmwareProgress() public delegate void Mcs::Usb::OnUpdateFirmwareProgress
(
             int )
10.2.2.11 OnUpdateFirmwareStatusChange() public delegate void Mcs::Usb::OnUpdateFirmware←
StatusChange (
            String^{\wedge} )
10.2.2.12 RoboStatusEventDelegate() public delegate void Mcs::Usb::RoboStatusEventDelegate (
             array< unsigned char >^{\wedge} buffer )
```

11 Class Documentation

11.1 CW2100_FunctionNet::AudioChannelsNet Struct Reference

Public Attributes

- W2100DacqGroupChannelEnumNet dacqgroup
- int channel
- · int amplification

11.1.1 Member Data Documentation

```
11.1.1.1 amplification int amplification
```

```
11.1.1.2 channel int channel
```

```
11.1.1.3 dacqgroup W2100DacqGroupChannelEnumNet dacqgroup
```

11.2 BatteryState Class Reference

Properties

- unsigned int Charge [get]
- unsigned int Voltage [get]
- System::String^ ChargeString [get]
- System::String^ ChargeRegionString [get]
- System::String^ VoltageString [get]

11.2.1 Property Documentation

```
11.2.1.1 Charge unsigned int Charge [get]
```

$\textbf{11.2.1.2} \quad \textbf{ChargeRegionString} \quad \texttt{System::} \quad \texttt{String}^{\wedge} \quad \texttt{ChargeRegionString} \quad \texttt{[get]}$


```
11.2.1.3 ChargeString System:: String^ ChargeString [get]
```

11.2.1.4 Voltage unsigned int Voltage [get]

11.2.1.5 VoltageString System:: String^ VoltageString [get]

11.3 BesselFilterBandPassNet Class Reference

Inheritance diagram for BesselFilterBandPassNet:

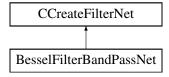

Public Member Functions

• BesselFilterBandPassNet (int numCoefSets, int order, double sampleRate, double centerFrequency, double widthFrequency, double scale)

Additional Inherited Members


11.3.1 Constructor & Destructor Documentation

11.3.1.1 BesselFilterBandPassNet() BesselFilterBandPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double centerFrequency,
double widthFrequency,
double scale )
```

11.4 BesselFilterHighPassNet Class Reference

Inheritance diagram for BesselFilterHighPassNet:

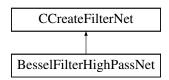

Public Member Functions

• BesselFilterHighPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.4.1 Constructor & Destructor Documentation

11.4.1.1 BesselFilterHighPassNet() BesselFilterHighPassNet (int numCoefSets,

```
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.5 BesselFilterLowPassNet Class Reference

Inheritance diagram for BesselFilterLowPassNet:

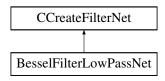

Public Member Functions

• BesselFilterLowPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.5.1 Constructor & Destructor Documentation

11.5.1.1 BesselFilterLowPassNet() BesselFilterLowPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.6 ButterworthFilterBandPassNet Class Reference

Inheritance diagram for ButterworthFilterBandPassNet:

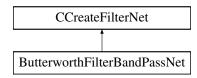

Public Member Functions

• ButterworthFilterBandPassNet (int numCoefSets, int order, double sampleRate, double centerFrequency, double widthFrequency, double scale)

Additional Inherited Members


11.6.1 Constructor & Destructor Documentation

11.6.1.1 ButterworthFilterBandPassNet() ButterworthFilterBandPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double centerFrequency,
double widthFrequency,
double scale )
```

11.7 ButterworthFilterHighPassNet Class Reference

Inheritance diagram for ButterworthFilterHighPassNet:

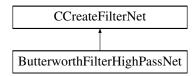

Public Member Functions

ButterworthFilterHighPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.7.1 Constructor & Destructor Documentation

11.7.1.1 ButterworthFilterHighPassNet() ButterworthFilterHighPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.8 ButterworthFilterLowPassNet Class Reference

Inheritance diagram for ButterworthFilterLowPassNet:

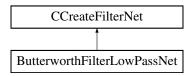

Public Member Functions

ButterworthFilterLowPassNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)

Additional Inherited Members


11.8.1 Constructor & Destructor Documentation

11.8.1.1 ButterworthFilterLowPassNet() ButterworthFilterLowPassNet (

```
int numCoefSets,
int order,
double sampleRate,
double cutoffFrequency,
double scale )
```

11.9 CChannelTestDeviceNet Class Reference

Inheritance diagram for CChannelTestDeviceNet:

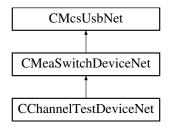

Public Member Functions

- CChannelTestDeviceNet ()
- ∼CChannelTestDeviceNet ()
- void SetWaveform (unsigned int Waveform)
- void SetAmplitude (unsigned int Amplitude)
- void SetFrequency (unsigned int Frequency)
- void SetAttenuation (unsigned int Attenuation)

Additional Inherited Members

11.9.1 Constructor & Destructor Documentation

```
11.9.1.1 CChannelTestDeviceNet() CChannelTestDeviceNet ( )
```

```
11.9.1.2 ~CChannelTestDeviceNet() ~CChannelTestDeviceNet ()
```


11.9.2 Member Function Documentation

```
11.9.2.1 SetAmplitude() void SetAmplitude (
unsigned int Amplitude )
```

```
11.9.2.3 SetFrequency() void SetFrequency (
          unsigned int Frequency )
```

11.10 CCMOSMea_FunctionNet Class Reference

Inheritance diagram for CCMOSMea_FunctionNet:

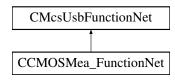

Public Member Functions

- CCMOSMea_FunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] cMOSMea_←
 FunctionPointerContainer)
- CCMOSMea_FunctionNet (CMcsUsbNet[^] mcsusb)
- void SetADCInputOffset (int32 t offset)
- int32_t GetADCInputOffset ()
- void SetSourceDrain (int32 t voltage)
- int32_t GetSourceDrain ()
- void SetSourceGate (int32_t voltage)
- int32_t GetSourceGate ()
- void SetSourceBulk (int32_t voltage)
- int32 t GetSourceBulk ()
- void SetGate (int32_t voltage)
- int32_t GetGate ()
- void SetBath (int32_t voltage)
- int32_t GetBath ()
- int32_t GetGNDI ()
- int32_t GetVDDI ()
- int32 t GetVDD3I ()
- void UpdateTransistorVoltages ()
- bool AreTransistorVoltagesSet ()
- void PowerChip (bool on)
- bool IsChipPowered ()
- enCMosMeaChipType DetectChipType ()
- void SetGateToVOP ()
- · void SetGateFloating ()
- bool IsGateFloating ()
- void VOPSTimerSetResetTimes (uint32_t ResetTime, uint32_t IntervalTime)
- void VOPSTimerSetResetTimes (uint32_t ResetTime, uint32_t IntervalTime, uint32_t HPFilterResetTime)
- void SetBathMode (CMOSMeaBathModeEnumNet Mode)
- CMOSMeaBathModeEnumNet GetBathMode ()
- void SetNeurochipMemoryData (uint16_t MemAddress, uint32_t MemData)
- void SetNeurochipMemoryData (uint16_t MemAddress, array< uint32_t >^ MemData)
- uint32 t GetNeurochipMemoryData (uint16 t MemAddress)
- $\bullet \ \, array < uint32_t > \ \, ^{\wedge} \, \, GetNeurochipMemoryData \, (uint16_t \, MemAddress, \, uint32_t \, ReqestLength)$
- uint32_t GetNeurochipMemorySize ()
- uint32_t GetMaxNumOfColumns (uint32_t Samplerate)
- void SetStimulusSites (List< int16_t >^ SwitchPosition)
- List< int16 t > ^ GetStimulusSites ()
- void ClearSTGOutput (uint32_t Channel)
- uint32_t GetNumberOfSupportedGroups ()
- uint32_t GetNumberOfSupportedGroups (uint32_t virtualDevice)

- DacqGroupChannelEnumNet GetGroupID (uint32_t Index)
- DacqGroupChannelEnumNet GetGroupID (uint32_t Index, uint32_t virtualDevice)
- uint32 t GetGroupNumberOfChannels (DacqGroupChannelEnumNet GroupID)
- uint32 t GetGroupNumberOfChannels (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- DacqMeaGroupTypeEnumNet GetGroupType (DacqGroupChannelEnumNet GroupID)
- void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >[^] EnabledChannelsBit←
 Map)
- void EnableChannelsInGroup (DacqGroupChannelEnumNet GroupID, List< bool >^ EnabledChannelsBit←
 Map, uint32 t virtualDevice)
- List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumNet GroupID)
- List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- uint32 t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID)
- uint32_t GetGroupResolutionPerDigit (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID)
- CMOSMeaValueUnitEnumNet GetGroupUnit (DacqGroupChannelEnumNet GroupID, uint32_t virtualDevice)
- int32_t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID)
- int32 t GetGroupDCOffset (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- int32_t GetGroupADCBits (DacqGroupChannelEnumNet GroupID)
- int32 t GetGroupADCBits (DacqGroupChannelEnumNet GroupID, uint32 t virtualDevice)
- uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel← Number)
- uint32_t GetGroupChannelBitmaskBySelect (DacqGroupChannelEnumNet GroupID, uint32_t Channel
 — Number, uint32_t virtualDevice)
- CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber)
- CMOSMeaInterfaceADCEnumNet GetGroupChannelBitmaskInterfaceADC (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMealFDigChannelEnumNet GetGroupChannelBitmaskIFDigChannels (uint32_t ChannelNumber)
- CMOSMealFDigChannelEnumNet GetGroupChannelBitmasklFDigChannels (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCBathCurrentEnumNet GetGroupChannelBitmaskHS1NCBathCurrent (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber)
- CMOSMeaHeadstage1NCCol2CurrentEnumNet GetGroupChannelBitmaskHS1NCCol2Current (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel
 — Number)
- CMOSMeaHeadstage1NChipTempEnumNet GetGroupChannelBitmaskHS1NChipTemp (uint32_t Channel
 — Number, uint32_t virtualDevice)
- CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t ChannelNumber)
- CMOSMeaSTG1DACSignalEnumNet GetGroupChannelBitmaskSTG1DACSignal (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32 t ChannelNumber)
- CMOSMeaHS1SidebandEnumNet GetGroupChannelBitmaskHS1Sidebands (uint32_t ChannelNumber, uint32_t virtualDevice)
- CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel
 — Number)
- CMOSMeaHS1TriggerStatusEnumNet GetGroupChannelBitmaskHS1TriggerStatus (uint32_t Channel↔ Number, uint32_t virtualDevice)

- CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_← t ChannelNumber)
- CMOSMeaPacketFrameContextGroupEnumNet GetGroupChannelBitmaskPacketFrameContext (uint32_← t ChannelNumber, uint32_t virtualDevice)

Additional Inherited Members

11.10.1 Constructor & Destructor Documentation

```
11.10.1.1 CCMOSMea_FunctionNet() [1/2] CCMOSMea_FunctionNet (
              CMcsUsbNet<sup>∧</sup> mcsusb,
              {\tt CMcsUsbFunctionPointerContainer}^{\land} \  \, {\tt cMOSMea\_FunctionPointerContainer} \  \, )
11.10.1.2 CCMOSMea_FunctionNet() [2/2] CCMOSMea_FunctionNet (
              CMcsUsbNet^ mcsusb )
11.10.2 Member Function Documentation
11.10.2.1 AreTransistorVoltagesSet() bool AreTransistorVoltagesSet ( )
11.10.2.2 ClearSTGOutput() void ClearSTGOutput (
              uint32_t Channel )
11.10.2.3 DetectChipType() enCMosMeaChipType DetectChipType ()
11.10.2.4 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (
              DacqGroupChannelEnumNet GroupID,
              List<br/>< bool >^{\land} EnabledChannelsBitMap )
```

```
11.10.2.5 EnableChannelsInGroup() [2/2] void EnableChannelsInGroup (
               DacqGroupChannelEnumNet GroupID,
               \label{eq:list_def} \mbox{List< bool} > ^{\wedge} \mbox{\it EnabledChannelsBitMap,}
               uint32_t virtualDevice )
11.10.2.6 GetADCInputOffset() int32_t GetADCInputOffset ()
11.10.2.7 GetBath() int32_t GetBath ()
11.10.2.8 GetBathMode() CMOSMeaBathModeEnumNet GetBathMode ( )
11.10.2.9 GetEnabledChannelsInGroup() [1/2] List<br/>bool> ^ GetEnabledChannelsInGroup (
               DacqGroupChannelEnumNet GroupID )
\textbf{11.10.2.10} \quad \textbf{GetEnabledChannelsInGroup()} \  \  \texttt{[2/2]} \quad \texttt{List} < \texttt{bool} > \  \  ^{\land} \  \  \texttt{GetEnabledChannelsInGroup} \  \  (
               DacqGroupChannelEnumNet GroupID,
               uint32_t virtualDevice )
11.10.2.11 GetGate() int32_t GetGate ()
11.10.2.12 GetGNDI() int32_t GetGNDI ( )
\textbf{11.10.2.13} \quad \textbf{GetGroupADCBits() [1/2]} \quad \texttt{int32\_t GetGroupADCBits ()}
               DacqGroupChannelEnumNet GroupID )
11.10.2.14 GetGroupADCBits() [2/2] int32_t GetGroupADCBits (
               DacqGroupChannelEnumNet GroupID,
               uint32_t virtualDevice )
```

```
11.10.2.15 GetGroupChannelBitmaskBySelect() [1/2] uint32_t GetGroupChannelBitmaskBySelect (
             DacqGroupChannelEnumNet GroupID,
             uint32_t ChannelNumber )
11.10.2.16 GetGroupChannelBitmaskBySelect() [2/2] uint32_t GetGroupChannelBitmaskBySelect (
             DacqGroupChannelEnumNet GroupID,
             uint32_t ChannelNumber,
             uint32_t virtualDevice )
11.10.2.17 GetGroupChannelBitmaskHS1NCBathCurrent() [1/2] CMOSMeaHeadstage1NCBathCurrentEnumNet
GetGroupChannelBitmaskHS1NCBathCurrent (
             uint32_t ChannelNumber )
11.10.2.18 GetGroupChannelBitmaskHS1NCBathCurrent() [2/2] CMOSMeaHeadstage1NCBathCurrentEnumNet
GetGroupChannelBitmaskHS1NCBathCurrent (
             uint32_t ChannelNumber,
             uint32_t virtualDevice )
11.10.2.19 GetGroupChannelBitmaskHS1NCCol2Current() [1/2] CMOSMeaHeadstage1NCCol2CurrentEnumNet
GetGroupChannelBitmaskHS1NCCol2Current (
             uint32_t ChannelNumber )
11.10.2.20 GetGroupChannelBitmaskHS1NCCol2Current() [2/2] CMOSMeaHeadstage1NCCol2CurrentEnumNet
GetGroupChannelBitmaskHS1NCCol2Current (
            uint32_t ChannelNumber,
             uint32_t virtualDevice )
11.10.2.21 GetGroupChannelBitmaskHS1NChipTemp() [1/2] CMOSMeaHeadstage1NChipTempEnumNet
GetGroupChannelBitmaskHS1NChipTemp (
            uint32_t ChannelNumber )
11.10.2.22 GetGroupChannelBitmaskHS1NChipTemp() [2/2] CMOSMeaHeadstage1NChipTempEnumNet
GetGroupChannelBitmaskHS1NChipTemp (
            uint32_t ChannelNumber,
             uint32_t virtualDevice )
```

```
11.10.2.23 GetGroupChannelBitmaskHS1Sidebands() [1/2] CMOSMeaHS1SidebandEnumNet GetGroup↔
ChannelBitmaskHS1Sidebands (
            uint32_t ChannelNumber )
11.10.2.24 GetGroupChannelBitmaskHS1Sidebands() [2/2] CMOSMeaHS1SidebandEnumNet GetGroup↔
ChannelBitmaskHS1Sidebands (
            uint32_t ChannelNumber,
            uint32_t virtualDevice )
11.10.2.25 GetGroupChannelBitmaskHS1TriggerStatus() [1/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
GroupChannelBitmaskHS1TriggerStatus (
            uint32_t ChannelNumber )
11.10.2.26 GetGroupChannelBitmaskHS1TriggerStatus() [2/2] CMOSMeaHS1TriggerStatusEnumNet Get↔
GroupChannelBitmaskHS1TriggerStatus (
             uint32_t ChannelNumber,
             uint32_t virtualDevice )
11.10.2.27 GetGroupChannelBitmasklFDigChannels() [1/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
ChannelBitmaskIFDigChannels (
            uint32_t ChannelNumber )
11.10.2.28 GetGroupChannelBitmaskIFDigChannels() [2/2] CMOSMeaIFDigChannelEnumNet GetGroup↔
ChannelBitmaskIFDigChannels (
            uint32_t ChannelNumber,
            uint32_t virtualDevice )
11.10.2.29 GetGroupChannelBitmaskInterfaceADC() [1/2] CMOSMeaInterfaceADCEnumNet GetGroup←
ChannelBitmaskInterfaceADC (
            uint32_t ChannelNumber )
11.10.2.30 GetGroupChannelBitmaskInterfaceADC() [2/2] CMOSMeaInterfaceADCEnumNet GetGroup↔
ChannelBitmaskInterfaceADC (
            uint32_t ChannelNumber,
             uint32_t virtualDevice )
```

```
11.10.2.31 GetGroupChannelBitmaskPacketFrameContext() [1/2] CMOSMeaPacketFrameContextGroupEnumNet
GetGroupChannelBitmaskPacketFrameContext (
                                            uint32_t ChannelNumber )
\textbf{11.10.2.32} \quad \textbf{GetGroupChannelBitmaskPacketFrameContext()} \quad \texttt{[2/2]} \quad \texttt{CMOSMeaPacketFrameContextGroupEnumNet}
{\tt GetGroupChannelBitmaskPacketFrameContext} \ \ (
                                             uint32_t ChannelNumber,
                                            uint32_t virtualDevice )
\textbf{11.10.2.33} \quad \textbf{GetGroupChannelBitmaskSTG1DACSignal()} \  \, \texttt{[1/2]} \quad \texttt{CMOSMeaSTG1DACSignalEnumNet} \quad \texttt{Get} \\ \leftarrow \quad \text{CMOSMeaSTG1DACSignalEnumNet} \\ \text{Comparison} \\ \text{Co
GroupChannelBitmaskSTG1DACSignal (
                                            uint32_t ChannelNumber )
11.10.2.34 GetGroupChannelBitmaskSTG1DACSignal() [2/2] CMOSMeaSTG1DACSignalEnumNet Get ←
GroupChannelBitmaskSTG1DACSignal (
                                             uint32_t ChannelNumber,
                                             uint32_t virtualDevice )
11.10.2.35 GetGroupDCOffset() [1/2] int32_t GetGroupDCOffset (
                                             DacqGroupChannelEnumNet GroupID )
11.10.2.36 GetGroupDCOffset() [2/2] int32_t GetGroupDCOffset (
                                             DacqGroupChannelEnumNet GroupID,
                                             uint32_t virtualDevice )
11.10.2.37 GetGroupID() [1/2] DacqGroupChannelEnumNet GetGroupID (
                                             uint32_t Index )
11.10.2.38 GetGroupID() [2/2] DacqGroupChannelEnumNet GetGroupID (
                                             uint32_t Index,
                                             uint32_t virtualDevice )
```

```
11.10.2.39 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumNet GroupID )
11.10.2.40 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.10.2.41 GetGroupResolutionPerDigit() [1/2] uint32_t GetGroupResolutionPerDigit (
             DacqGroupChannelEnumNet GroupID )
11.10.2.42 GetGroupResolutionPerDigit() [2/2] uint32_t GetGroupResolutionPerDigit (
             DacqGroupChannelEnumNet GroupID,
            uint32_t virtualDevice )
11.10.2.43 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumNet GroupID )
11.10.2.44 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.10.2.45 GetGroupType() [1/2] DacqMeaGroupTypeEnumNet GetGroupType (
             DacqGroupChannelEnumNet GroupID )
11.10.2.46 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
             DacqGroupChannelEnumNet GroupID,
             uint32_t virtualDevice )
11.10.2.47 GetGroupUnit() [1/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
             DacqGroupChannelEnumNet GroupID )
```

```
11.10.2.48 GetGroupUnit() [2/2] CMOSMeaValueUnitEnumNet GetGroupUnit (
               DacqGroupChannelEnumNet GroupID,
               uint32_t virtualDevice )
11.10.2.49 GetMaxNumOfColumns() uint32_t GetMaxNumOfColumns (
               uint32_t Samplerate )
11.10.2.50 GetNeurochipMemoryData() [1/2] uint32_t GetNeurochipMemoryData (
               uint16_t MemAddress )
11.10.2.51 GetNeurochipMemoryData() [2/2] array<uint32_t> ^ GetNeurochipMemoryData (
               uint16_t MemAddress,
               uint32_t ReqestLength )
11.10.2.52 GetNeurochipMemorySize() uint32_t GetNeurochipMemorySize ( )
\textbf{11.10.2.53} \quad \textbf{GetNumberOfSupportedGroups()} \  \, \textbf{[1/2]} \quad \text{uint} \  \, \textbf{32\_t} \  \, \textbf{GetNumberOfSupportedGroups} \  \, \textbf{()} \\
\textbf{11.10.2.54} \quad \textbf{GetNumberOfSupportedGroups()} \  \  \texttt{[2/2]} \quad \texttt{uint32\_t} \  \  \texttt{GetNumberOfSupportedGroups} \  \  \textbf{(}
               uint32_t virtualDevice )
11.10.2.55 GetSourceBulk() int32_t GetSourceBulk ( )
11.10.2.56 GetSourceDrain() int32_t GetSourceDrain ( )
11.10.2.57 GetSourceGate() int32_t GetSourceGate ( )
```

```
11.10.2.58 GetStimulusSites() List<intl6_t> ^{\land} GetStimulusSites ( )
11.10.2.59 GetVDD3I() int32_t GetVDD3I ( )
11.10.2.60 GetVDDI() int32_t GetVDDI ( )
11.10.2.61 IsChipPowered() bool IsChipPowered ( )
11.10.2.62 IsGateFloating() bool IsGateFloating ( )
11.10.2.63 PowerChip() void PowerChip (
            bool on )
11.10.2.64 SetADCInputOffset() void SetADCInputOffset (
             int32_t offset )
11.10.2.65 SetBath() void SetBath (
            int32_t voltage )
11.10.2.66 SetBathMode() void SetBathMode (
            CMOSMeaBathModeEnumNet Mode )
11.10.2.67 SetGate() void SetGate (
            int32_t voltage )
```


```
11.10.2.68 SetGateFloating() void SetGateFloating ( )
11.10.2.69 SetGateToVOP() void SetGateToVOP ()
11.10.2.70 SetNeurochipMemoryData() [1/2] void SetNeurochipMemoryData (
             uint16_t MemAddress,
             array< uint32_t >^{\land} MemData )
11.10.2.71 SetNeurochipMemoryData() [2/2] void SetNeurochipMemoryData (
             uint16_t MemAddress,
             uint32_t MemData )
11.10.2.72 SetSourceBulk() void SetSourceBulk (
             int32_t voltage )
11.10.2.73 SetSourceDrain() void SetSourceDrain (
             int32_t voltage )
11.10.2.74 SetSourceGate() void SetSourceGate (
             int32_t voltage )
11.10.2.75 SetStimulusSites() void SetStimulusSites (
             List< int16_t >^{\land} SwitchPosition )
11.10.2.76 UpdateTransistorVoltages() void UpdateTransistorVoltages ( )
11.10.2.77 VOPSTimerSetResetTimes() [1/2] void VOPSTimerSetResetTimes (
             uint32_t ResetTime,
             uint32_t IntervalTime )
```

11.10.2.78 VOPSTimerSetResetTimes() [2/2] void VOPSTimerSetResetTimes (

```
uint32_t ResetTime,
uint32_t IntervalTime,
uint32_t HPFilterResetTime )
```

11.11 CCMOSMeaDeviceNet Class Reference

Inheritance diagram for CCMOSMeaDeviceNet:

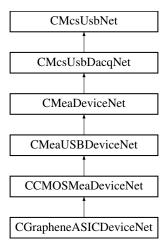

Classes

class CRegionOfInterestRect

Public Member Functions

- CCMOSMeaDeviceNet (void)
- ∼CCMOSMeaDeviceNet ()
- virtual void SetBaseSamplerate (int BaseSamplerate)
- int GetBaseSamplerate ()
- virtual array< int > ^ GetAvailableBaseSamplerates ()
- int GetMaxReadableColumns ()
- virtual void SetRegionOfInterests (System::Collections::Generic::Dictionary< int, CRegionOfInterestRect^>^
 rois)
- void UpdateChannelBlock (int queuesize, int threshold, int channels_in_block)
- System::Collections::Generic::Dictionary< int, array< array< int16_t >^>^> ^ GetCMOSDataDictionary
 (int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ GetChannelDataUI16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames ret)
- System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ GetChannelDataI16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< uint32_t>^> ^ GetChannelDataUI32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)
- System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ GetChannelDatal32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Properties

```
• CCMOSMea FunctionNet^ CMosMea [get]
```

```
• CStimulusFunctionNet^ Stimulus [get]
```

Additional Inherited Members

11.11.1 Constructor & Destructor Documentation

```
11.11.2.5 GetChannelDataUl16() System::Collections::Generic::Dictionary<int, array<uint16_←
t>^> ^{\wedge} GetChannelDataUI16 (
             DacqGroupChannelEnumNet group,
             int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.11.2.6 GetChannelDataUl32() System::Collections::Generic::Dictionary<int, array<uint32_←
t>^{\wedge}> ^ GetChannelDataUI32 (
             DacqGroupChannelEnumNet group,
             int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.11.2.7 GetCMOSDataDictionary() System::Collections::Generic::Dictionary<int, array<array<int16↔
int frames,
             [System::Runtime::InteropServices::Out] int % frames_ret )
11.11.2.8 GetMaxReadableColumns() int GetMaxReadableColumns ()
11.11.2.9 SetBaseSamplerate() virtual void SetBaseSamplerate (
             int BaseSamplerate ) [virtual]
Reimplemented in CGrapheneASICDeviceNet.
11.11.2.10 SetRegionOfInterests() virtual void SetRegionOfInterests (
             {\tt System::Collections::Generic::Dictionary<~int,~CRegionOfInterestRect}^{\wedge}{\gt}^{\wedge}~\textit{rois}~)
[virtual]
11.11.2.11 UpdateChannelBlock() void UpdateChannelBlock (
             int queuesize,
             int threshold,
             int channels_in_block )
```


11.11.3 Property Documentation

11.11.3.1 CMosMea CCMOSMea_FunctionNet^ CMosMea [get]

11.11.3.2 Stimulus CStimulusFunctionNet^ Stimulus [get]

11.12 CCreateFilterNet Class Reference

Inheritance diagram for CCreateFilterNet:

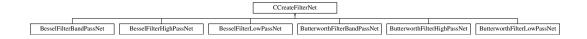

Public Member Functions

- · CCreateFilterNet (int numCoefSets, int order, double sampleRate, double cutoffFrequency, double scale)
- ∼CCreateFilterNet ()
- CFilterCoefficientsNet ^ GetBiQuad (int index)
- array< CFilterCoefficientsNet[^]> [^] GetBiQuads ()
- CCreateFilter * GetCpp ()

Static Public Member Functions

- static int FindFilter (array < CFilterCoefficientsNet[^]> coef, array < CCreateFilterNet[^]> param)
- static int FindFilter (array< array< uint64_t >^> coef, array< CCreateFilterNet^> param, CFilterCoefficientsNet::s_FilterAttributesNet^ FiltAttr, bool DoMCSLegacyCompare)

Protected Member Functions

• CCreateFilterNet (int numCoefSets, CCreateFilter *pCreateFilter)

Properties

- int NumCoefSets [get]
- int Order [get]
- double SampleRate [get]
- double CutoffOrCenterFrequency [get]
- double WidthFrequency [get]
- double Scale [get]

11.12.1 Constructor & Destructor Documentation

```
11.12.1.1 CCreateFilterNet() [1/2] CCreateFilterNet (
             int numCoefSets,
             int order,
             double sampleRate,
             double cutoffFrequency,
             double scale )
11.12.1.2 ~CCreateFilterNet() ~CCreateFilterNet ()
11.12.1.3 CCreateFilterNet() [2/2] CCreateFilterNet (
             int numCoefSets,
             CCreateFilter * pCreateFilter ) [protected]
11.12.2 Member Function Documentation
11.12.2.1 FindFilter() [1/2] static int FindFilter (
             array< array< uint64_t >^{\wedge}>^{\wedge} coef,
             array< CCreateFilterNet^>^ param,
             CFilterCoefficientsNet::s_FilterAttributesNet^ FiltAttr,
             bool DoMCSLegacyCompare ) [static]
11.12.2.2 FindFilter() [2/2] static int FindFilter (
             \textbf{11.12.2.3} \quad \textbf{GetBiQuad()} \quad \texttt{CFilterCoefficientsNet} \; \wedge \; \texttt{GetBiQuad} \; \; \textbf{(}
             int index )
11.12.2.4 GetBiQuads() array<CFilterCoefficientsNet^> ^ GetBiQuads ()
11.12.2.5 GetCpp() CCreateFilter* GetCpp ()
```

11.12.3 Property Documentation

```
11.12.3.1 CutoffOrCenterFrequency double CutoffOrCenterFrequency [get]
```

```
11.12.3.2 NumCoefSets int NumCoefSets [get]
```

```
11.12.3.3 Order int Order [get]
```


11.12.3.4 SampleRate double SampleRate [get]

11.12.3.5 Scale double Scale [get]

11.12.3.6 WidthFrequency double WidthFrequency [get]

11.13 CDacCalibrationFunctionNet Class Reference

Inheritance diagram for CDacCalibrationFunctionNet:

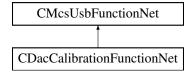

Public Member Functions

CDacCalibrationFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pDac
 — CalibrationFunctionPointerContainer)

Initializes a new instance of the CDacCalibrationFunctionNet class.

- CDacCalibrationFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CDacCalibrationFunctionNet ()
- !CDacCalibrationFunctionNet ()
- void SetDacOffset (uint16_t dacChannel, int32_t offset)

Sets the offset of a DAC channel.

int32_t GetDacOffset (uint16_t dacChannel)

Gets the offset of a DAC channel.

void BurnDacOffset (uint16_t dacChannel)

Writes the offset of a DAC channel to permanent memory.

Additional Inherited Members

11.13.1 Detailed Description

11.13.2 Constructor & Destructor Documentation

```
11.13.2.1 CDacCalibrationFunctionNet() [1/2] CDacCalibrationFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pDacCalibrationFunctionPointerContainer)
```

Initializes a new instance of the CDacCalibrationFunctionNet class.

```
11.13.2.2 CDacCalibrationFunctionNet() [2/2] CDacCalibrationFunctionNet ( CMcsUsbNet^ mcsusb )
```

```
11.13.2.3 ~CDacCalibrationFunctionNet() virtual ~CDacCalibrationFunctionNet ( ) [virtual]
```

```
11.13.2.4 "!CDacCalibrationFunctionNet() !CDacCalibrationFunctionNet ( )
```

11.13.3 Member Function Documentation

```
11.13.3.1 BurnDacOffset() void BurnDacOffset ( uint16_t dacChannel )
```

Writes the offset of a DAC channel to permanent memory.

Parameters

```
dacChannel The DAC channel number.
```

Gets the offset of a DAC channel.

Parameters

Returns

The offset in digits.


Sets the offset of a DAC channel.

Parameters

dacChannel	The DAC channel number.
offset	The offset in digits.

11.14 CDacqGroupChannelGenericSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelGenericSelectionNet:

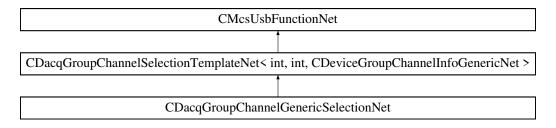

Public Member Functions

CDacqGroupChannelGenericSelectionNet (CMcsUsbNet[^] mcsusb)


Additional Inherited Members

11.14.1 Constructor & Destructor Documentation

11.14.1.1 CDacqGroupChannelGenericSelectionNet() CDacqGroupChannelGenericSelectionNet (CMcsUsbNet^ mcsusb)

11.15 CDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CDacqGroupChannelSelectionNet:

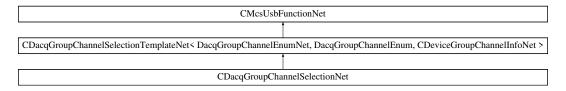

Public Member Functions

CDacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members

11.15.1 Constructor & Destructor Documentation

```
11.15.1.1 CDacqGroupChannelSelectionNet() CDacqGroupChannelSelectionNet (
CMcsUsbNet^ mcsusb )
```

11.16 CDacqGroupChannelSelectionTemplateNet < DacqGroupChannelEnumTemplateNet, DacqGroupChannelEnumTemplate, CDeviceGroupChannelInfoTemplateNet > Class Template Reference

 $Inheritance \ diagram \ for \ CDacqGroupChannelSelectionTemplateNet < \ DacqGroupChannelEnumTemplateNet, \\ DacqGroupChannelEnumTemplate, \ CDeviceGroupChannelInfoTemplateNet >: \\$

```
CMcsUsbFunctionNet

CDacqGroupChannelSelectionTemplateNet < DacqGroupChannelEnumTemplateNet > DacqGroupChann
```

Public Member Functions

- CDacqGroupChannelSelectionTemplateNet (CMcsUsbNet[^] mcsusb)
- uint32 t GetNumberOfSupportedGroups ()
- uint32_t GetNumberOfSupportedGroups (uint32_t virtualDevice)
- DacqGroupChannelEnumTemplateNet GetGroupID (uint32_t Index)
- DacqGroupChannelEnumTemplateNet GetGroupID (uint32_t Index, uint32_t virtualDevice)
- uint32_t GetGroupNumberOfChannels (DacqGroupChannelEnumTemplateNet GroupID)
- DacgMeaGroupTypeEnumNet GetGroupType (DacgGroupChannelEnumTemplateNet GroupID)
- DacqMeaGroupTypeEnumNet GetGroupType (DacqGroupChannelEnumTemplateNet GroupID, uint32_

 t virtualDevice)

- void EnableChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, List< bool >[^] Enabled←
 ChannelsBitMap)
- void EnableChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, List< bool >^ Enabled←
 ChannelsBitMap, uint32_t virtualDevice)
- List< bool > ^ GetEnabledChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID)
- List< bool > [^] GetEnabledChannelsInGroup (DacqGroupChannelEnumTemplateNet GroupID, uint32_← t virtualDevice)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumTemplateNet GroupID)
- SampleSizeNet GetGroupSampleSize (DacqGroupChannelEnumTemplateNet GroupID, uint32_t virtual → Device)
- List< CDeviceGroupChannelInfoTemplateNet[^]> [^] GetDeviceGroupChannelInfos ()
- List< CDeviceGroupChannelInfoTemplateNet^> ^ GetDeviceGroupChannelInfos (uint32_t virtualDevice)

Additional Inherited Members

11.16.1 Constructor & Destructor Documentation

```
11.16.1.1 CDacqGroupChannelSelectionTemplateNet() CDacqGroupChannelSelectionTemplateNet (
CMcsUsbNet^ mcsusb )
```

11.16.2 Member Function Documentation

```
11.16.2.1 EnableChannelsInGroup() [1/2] void EnableChannelsInGroup (

DacqGroupChannelEnumTemplateNet GroupID,

List< bool >^ EnabledChannelsBitMap )
```

```
11.16.2.2 EnableChannelsInGroup() [2/2] void EnableChannelsInGroup (
DacqGroupChannelEnumTemplateNet GroupID,
```

```
List< bool > ^ EnabledChannelsBitMap, uint32_t virtualDevice)
```

```
11.16.2.5 GetEnabledChannelsInGroup() [1/2] List<bool> ^ GetEnabledChannelsInGroup (
             DacqGroupChannelEnumTemplateNet GroupID )
11.16.2.6 GetEnabledChannelsInGroup() [2/2] List<br/>bool> ^ GetEnabledChannelsInGroup (
             DacqGroupChannelEnumTemplateNet GroupID,
             uint32_t virtualDevice )
11.16.2.7 GetGroupID() [1/2] DacqGroupChannelEnumTemplateNet GetGroupID (
             uint32_t Index )
11.16.2.8 GetGroupID() [2/2] DacqGroupChannelEnumTemplateNet GetGroupID (
             uint32_t Index,
             uint32_t virtualDevice )
11.16.2.9 GetGroupNumberOfChannels() [1/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumTemplateNet GroupID )
11.16.2.10 GetGroupNumberOfChannels() [2/2] uint32_t GetGroupNumberOfChannels (
             DacqGroupChannelEnumTemplateNet GroupID,
             uint32_t virtualDevice )
11.16.2.11 GetGroupSampleSize() [1/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumTemplateNet GroupID )
11.16.2.12 GetGroupSampleSize() [2/2] SampleSizeNet GetGroupSampleSize (
             DacqGroupChannelEnumTemplateNet GroupID,
             uint32_t virtualDevice )
11.16.2.13 GetGroupType() [1/2] DacqMeaGroupTypeEnumNet GetGroupType (
             DacqGroupChannelEnumTemplateNet GroupID )
```


```
11.16.2.14 GetGroupType() [2/2] DacqMeaGroupTypeEnumNet GetGroupType (
DacqGroupChannelEnumTemplateNet GroupID,
uint32_t virtualDevice )

11.16.2.15 GetNumberOfSupportedGroups() [1/2] uint32_t GetNumberOfSupportedGroups ( )

11.16.2.16 GetNumberOfSupportedGroups() [2/2] uint32_t GetNumberOfSupportedGroups ( uint32_t virtualDevice )
```

11.17 CDeviceGroupChannelInfoGenericNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoGenericNet:

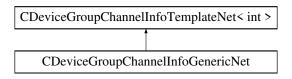

Public Member Functions


• CDeviceGroupChannelInfoGenericNet (int id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.17.1 Constructor & Destructor Documentation

11.18 CDeviceGroupChannelInfoMEA2100_256Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoMEA2100_256Net:

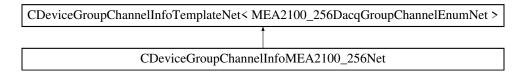

Public Member Functions


 CDeviceGroupChannelInfoMEA2100_256Net (MEA2100_256DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.18.1 Constructor & Destructor Documentation

11.19 CDeviceGroupChannelInfoNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoNet:

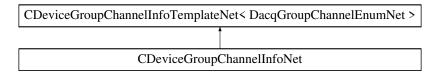

Public Member Functions


CDeviceGroupChannelInfoNet (DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.19.1 Constructor & Destructor Documentation

11.20 CDeviceGroupChannelInfoSCUNet Class Reference

Inheritance diagram for CDeviceGroupChannelInfoSCUNet:

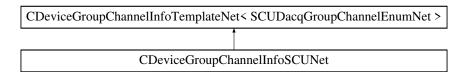

Public Member Functions

CDeviceGroupChannelInfoSCUNet (SCUDacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.20.1 Constructor & Destructor Documentation

${\bf 11.21} \quad {\bf CDevice Group Channel Info Template Net} < {\bf Dacq Group Channel Enum Template Net}$

> Class Template Reference

Public Member Functions

CDeviceGroupChannelInfoTemplateNet (DacqGroupChannelEnumTemplateNet id, int channels, DacqMeaGroupTypeEnumNet type)

Public Attributes

- DacqGroupChannelEnumTemplateNet GroupID
- int NumberOfChannels
- DacqMeaGroupTypeEnumNet GroupType

11.21.1 Constructor & Destructor Documentation

11.21.1.1 CDeviceGroupChannelInfoTemplateNet() CDeviceGroupChannelInfoTemplateNet (

```
DacqGroupChannelEnumTemplateNet id,
int channels,
DacqMeaGroupTypeEnumNet type )
```

11.21.2 Member Data Documentation


11.21.2.1 GroupID DacqGroupChannelEnumTemplateNet GroupID

11.21.2.2 GroupType DacqMeaGroupTypeEnumNet GroupType

11.21.2.3 NumberOfChannels int NumberOfChannels

11.22 CDeviceGroupChannelInfoW2100Net Class Reference

Inheritance diagram for CDeviceGroupChannelInfoW2100Net:

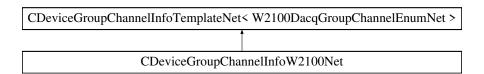

Public Member Functions

• CDeviceGroupChannelInfoW2100Net (W2100DacqGroupChannelEnumNet id, int channels, DacqMeaGroupTypeEnumNet type)

Additional Inherited Members

11.22.1 Constructor & Destructor Documentation


11.22.1.1 CDeviceGroupChannelInfoW2100Net() CDeviceGroupChannelInfoW2100Net (

```
W2100DacqGroupChannelEnumNet id,
int channels,
DacqMeaGroupTypeEnumNet type )
```

11.23 CDigOutStimulatorFunctionNet Class Reference

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

Inheritance diagram for CDigOutStimulatorFunctionNet:

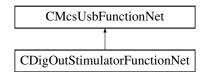

Public Member Functions

CDigOutStimulatorFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pDigOut
 — StimulatorFunctionPointerContainer)

Initializes a new instance of the CDigOutStimulatorFunctionNet class.

- CDigOutStimulatorFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CDigOutStimulatorFunctionNet ()
- !CDigOutStimulatorFunctionNet ()
- void ClearChannel (int32 t NrChannel)

clear stimulation pattern

CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareChannelData (array< int32_t >^ Amplitude, array< uint64_t >^ Duration)

prepares the channel data for the device and unrolles the data for the GUI

void SendChannelData (int32_t NrChannel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled)

send or append stimulation pattern

• int32 t GetNumberOfChannels ()

get the number of channels available on the device

• void SetGlobalRepeat (int32_t NrChannel, bool value)

set repeat whole stimulation pattern

bool GetGlobalRepeat (int32_t NrChannel)

get repeat whole stimulation pattern

• void SetStartTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)

sets start condition of digital out stimulator

DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope (int32_t NrChannel)

queries start condition of digital out stimulator

• void SetStopTriggerSlope (int32_t NrChannel, DigitalStimulatorTriggerSlopeEnumNet Condition)

sets stop condition of digital out stimulator

DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope (int32_t NrChannel)

queries stop condition of digital out stimulator

Additional Inherited Members

11.23.1 Detailed Description

CDigOutStimulatorFunctionNet is the class of the DigOut stimulator function class.

11.23.2 Constructor & Destructor Documentation

```
11.23.2.1 CDigOutStimulatorFunctionNet() [1/2] CDigOutStimulatorFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pDigOutStimulatorFunctionPointerContainer)
```

 $\label{lem:linear_continuous_co$

```
11.23.2.2 CDigOutStimulatorFunctionNet() [2/2] CDigOutStimulatorFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
\textbf{11.23.2.3} \quad \sim \textbf{CDigOutStimulatorFunctionNet()} \quad \text{virtual} \quad \sim \textbf{CDigOutStimulatorFunctionNet ()} \quad \text{[virtual]}
```

```
11.23.2.4 "!CDigOutStimulatorFunctionNet() !CDigOutStimulatorFunctionNet ( )
```

11.23.3 Member Function Documentation

```
11.23.3.1 ClearChannel() void ClearChannel (
    int32_t NrChannel )
```

clear stimulation pattern

Parameters

NrChannel the channel to clear

get repeat whole stimulation pattern

Parameters

NrChannel channel number

Returns

current value

```
11.23.3.3 GetNumberOfChannels() int32_t GetNumberOfChannels ( )
```

get the number of channels available on the device

Returns

the number of channels

```
11.23.3.4 GetStartTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetStartTriggerSlope ( int32_t NrChannel )
```

queries start condition of digital out stimulator

Parameters

NrChannel channel number

Returns

start condition (rising or falling edge)

```
11.23.3.5 GetStopTriggerSlope() DigitalStimulatorTriggerSlopeEnumNet GetStopTriggerSlope ( int32_t NrChannel )
```

queries stop condition of digital out stimulator

Parameters

```
NrChannel channel number
```

Returns

stop condition (rising or falling edge)

11.23.3.6 PrepareChannelData() CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ Prepare← ChannelData (

```
array< int32_t >^{\land} Amplitude, array< uint64_t >^{\land} Duration )
```

prepares the channel data for the device and unrolles the data for the GUI

Parameters

Amplitude	array of amplitudes
Duration	array of durations

Returns

send or append stimulation pattern

Parameters

NrChannel	the channel to send data to
device_data_and_unrolled	internal, use value obtained from PrepareChannelData

set repeat whole stimulation pattern

Parameters

NrChannel	channel number
value	new value

sets start condition of digital out stimulator

NrChannel	channel number
Condition	start condition (rising or falling edge)

sets stop condition of digital out stimulator


Parameters

NrChannel	channel number
Condition	stop condition (rising or falling edge)

11.24 CEncapsulatorDeviceNet Class Reference

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

Inheritance diagram for CEncapsulatorDeviceNet:

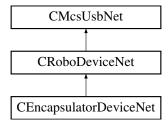

Public Member Functions

- CEncapsulatorDeviceNet (void)
- CRoboFluidDeviceNet ^ GetRoboFluidDevice ()

Additional Inherited Members

11.24.1 Detailed Description

CEncapsulatorDeviceNet is the to control the MCS HiClamp device

11.24.2 Constructor & Destructor Documentation

11.24.2.1 CEncapsulatorDeviceNet() CEncapsulatorDeviceNet (void)


11.24.3 Member Function Documentation

11.24.3.1 GetRoboFluidDevice() CRoboFluidDeviceNet $^{\land}$ GetRoboFluidDevice ()

11.25 CExternDTesterDeviceNet Class Reference

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

Inheritance diagram for CExternDTesterDeviceNet:

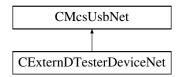

Public Member Functions

CExternDTesterDeviceNet ()

Initializes a new instance of the CExternDTesterDeviceNet class.

- virtual ~CExternDTesterDeviceNet ()
- !CExternDTesterDeviceNet ()
- array< uint8_t > ^ Read (int configString_Length)

Reads the config string from the device.

String ^ Read2 ()

Reads the config string from the device.

void Write (array< uint8_t >^ configString)

Reads the config string from the device.

void Write2 (String[^] configString)

Reads the config string from the device.

Additional Inherited Members

11.25.1 Detailed Description

CExternDTesterDeviceNet is the class to access the ExternD Tester (Handheld Device Tester D)

11.25.2 Constructor & Destructor Documentation

```
11.25.2.1 CExternDTesterDeviceNet() CExternDTesterDeviceNet ()
```

Initializes a new instance of the CExternDTesterDeviceNet class.

```
11.25.2.2 ~CExternDTesterDeviceNet() virtual ~CExternDTesterDeviceNet () [virtual]
```

```
11.25.2.3 "!CExternDTesterDeviceNet() !CExternDTesterDeviceNet ()
```

11.25.3 Member Function Documentation

Reads the config string from the device.

Parameters

Returns

The config string.

```
11.25.3.2 Read2() String ^{\wedge} Read2 ( )
```

Reads the config string from the device.

Returns

The config string.

```
11.25.3.3 Write() void Write (
array < uint8_t >^{\land} configString)
```

Reads the config string from the device.

configString	The config string.
--------------	--------------------

```
11.25.3.4 Write2() void Write2 (
String^ configString)
```

Reads the config string from the device.

Parameters

configString The config string.	configString	The config string.
---------------------------------	--------------	--------------------

11.26 CFilterCoefficientsNet Class Reference

Classes

struct s_FilterAttributesNet

Public Member Functions

- CFilterCoefficientsNet ()
- CFilterCoefficientsNet (double b0, double b1, double b2, double a1, double a2)
- CFilterCoefficientsNet (double b0, double b1, double a1)
- CFilterCoefficientsNet (array< double >[^] b, array< double >[^] a)
- ∼CFilterCoefficientsNet ()
- bool IsEqual (array< uint64_t >^ coefficients, s_FilterAttributesNet^ FiltAttr)
- bool IsEqual (array< uint64_t >^ coefficients, s_FilterAttributesNet^ FiltAttr, bool DoMCSLegacyCompare)
- uint64_t GetUintB (int index, s_FilterAttributesNet[^] FiltAttr)
- uint64_t GetUintA (int index, s_FilterAttributesNet[^] FiltAttr)

Properties

```
array< double > A [get]
array< double > B [get]
```

11.26.1 Constructor & Destructor Documentation

11.26.1.1 CFilterCoefficientsNet() [1/4] CFilterCoefficientsNet ()

```
11.26.1.2 CFilterCoefficientsNet() [2/4] CFilterCoefficientsNet (
              double b0,
              double b1,
              double b2,
              double a1,
              double a2 )
11.26.1.3 CFilterCoefficientsNet() [3/4] CFilterCoefficientsNet (
              double b0,
              double b1,
              double a1 )
11.26.1.4 CFilterCoefficientsNet() [4/4] CFilterCoefficientsNet (
              array< double >^{\wedge} b,
              array< double >^{\wedge} a)
11.26.1.5 ~CFilterCoefficientsNet() ~CFilterCoefficientsNet ()
11.26.2 Member Function Documentation
11.26.2.1 GetUintA() uint64_t GetUintA (
              int index,
              s_FilterAttributesNet^ FiltAttr )
\textbf{11.26.2.2} \quad \textbf{GetUintB()} \quad \texttt{uint64\_t GetUintB} \ \ (
              int index,
              s_FilterAttributesNet^ FiltAttr )
11.26.2.3 | IsEqual() [1/2] | bool IsEqual (
              array< uint64_t >^{\land} coefficients,
              s_FilterAttributesNet^ FiltAttr )
```


11.26.3 Property Documentation

```
11.26.3.1 A array< double>^ A [get]
```

11.26.3.2 B array< double> $^{\land}$ B [get]

11.27 CFilterConfigurationNet Class Reference

Inheritance diagram for CFilterConfigurationNet:

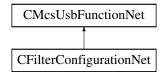

Public Member Functions

- CFilterConfigurationNet (CMcsUsbNet[^] mcsusb)
- void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet[^] Coefficients, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameter (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, CFilterCoefficientsNet[^] CoefficientsSet1, CFilterCoefficientsNet[^] CoefficientsSet2, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)
- void EraseFilterParameterPermanent (DacqGroupChannelEnumNet GroupID, uint32 t FilterNumber)
- · void SetHighpassFilterEnable (bool enable)
- bool GetHighpassFilterEnable ()
- void ResetHighpassFilter ()
- uint32_t GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber, FilterAttributeEnumNet index)
- CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilterAttributes (DacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)

Static Public Member Functions

static CFilterCoefficientsNet::s FilterAttributesNet ^ GetDefaultFilterAttributes ()

Additional Inherited Members

```
11.27.1 Constructor & Destructor Documentation
```

```
11.27.1.1 CFilterConfigurationNet() CFilterConfigurationNet (
             CMcsUsbNet^ mcsusb )
11.27.2 Member Function Documentation
11.27.2.1 EraseFilterParameterPermanent() void EraseFilterParameterPermanent (
             DacqGroupChannelEnumNet GroupID,
             uint32\_t FilterNumber )
11.27.2.2 GetDefaultFilterAttributes() static CFilterCoefficientsNet::s_FilterAttributesNet ^
GetDefaultFilterAttributes ( ) [static]
11.27.2.3 GetFilterAttributes() [1/2] CFilterCoefficientsNet::s_FilterAttributesNet ^ GetFilter←
Attributes (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber )
11.27.2.4 GetFilterAttributes() [2/2] uint32_t GetFilterAttributes (
             DacqGroupChannelEnumNet GroupID,
             uint32_t FilterNumber,
             FilterAttributeEnumNet index )
11.27.2.5 GetHighpassFilterEnable() bool GetHighpassFilterEnable ( )
```

11.27.2.6 ResetHighpassFilter() void ResetHighpassFilter ()

11.27.2.7 SetFilterParameter() [1/2] void SetFilterParameter (

```
DacqGroupChannelEnumNet GroupID,

uint32_t FilterNumber,

CFilterCoefficientsNet^ Coefficients,

CFilterPropertyNet^ FilterProp )

11.27.2.8 SetFilterParameter() [2/2] void SetFilterParameter (

DacqGroupChannelEnumNet GroupID,

uint32_t FilterNumber,

CFilterCoefficientsNet^ CoefficientsSet1,

CFilterCoefficientsNet^ CoefficientsSet2,

CFilterPropertyNet^ FilterProp )

11.27.2.9 SetFilterParameterPermanent() void SetFilterParameterPermanent (


DacqGroupChannelEnumNet GroupID,

uint32_t FilterNumber )

11.27.2.10 SetHighpassFilterEnable() void SetHighpassFilterEnable (
bool enable )
```

11.28 CFilterConfigurationRegisterNet Class Reference

Inheritance diagram for CFilterConfigurationRegisterNet:

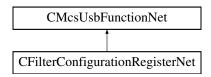

Public Member Functions

- CFilterConfigurationRegisterNet (CMcsUsbNet[^] mcsusb)
- void SetFilterParameter (uint32_t FilterCoefRegBase, CFilterCoefficientsNet[^] Coefficients, uint32_t Filter←
 InfoRegBase, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameter (uint32_t FilterCoefSet1RegBase, CFilterCoefficientsNet[^] CoefficientsSet1, uint32_t FilterCoefSet2RegBase, CFilterCoefficientsNet[^] CoefficientsSet2, uint32_t FilterInfoRegBase, CFilterPropertyNet[^] FilterProp)
- void SetFilterParameterPermanent (uint32_t FilterCoefRegBase, uint32_t FilterCoefDmaReg, uint32_← t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBize)
- void SetFilterParameterPermanent (uint32_t FilterCoefSet1RegBase, uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2RegBase, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoRegBase, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefDmaReg, uint32_t FilterInfoDmaReg, uint32_t EEP-ROMBase, uint32_t EEPROMSize)
- void EraseFilterParameterPermanent (uint32_t FilterCoefSet1DmaReg, uint32_t FilterCoefSet2DmaReg, uint32_t FilterInfoDmaReg, uint32_t EEPROMBase, uint32_t EEPROMSize)

Additional Inherited Members

11.28.1 Constructor & Destructor Documentation

```
11.28.1.1 CFilterConfigurationRegisterNet() CFilterConfigurationRegisterNet ( CMcsUsbNet^ mcsusb )
```

11.28.2 Member Function Documentation

11.28.2.2 EraseFilterParameterPermanent() [2/2] void EraseFilterParameterPermanent (

```
uint32_t FilterCoefSet1DmaReg,
uint32_t FilterCoefSet2DmaReg,
uint32_t FilterInfoDmaReg,
uint32_t EEPROMBase,
uint32_t EEPROMSize)
```

11.28.2.3 SetFilterParameter() [1/2] void SetFilterParameter (

```
uint32_t FilterCoefRegBase,
CFilterCoefficientsNet^ Coefficients,
uint32_t FilterInfoRegBase,
CFilterPropertyNet^ FilterProp )
```

11.28.2.4 SetFilterParameter() [2/2] void SetFilterParameter (

```
uint32_t FilterCoefSet1RegBase,
CFilterCoefficientsNet^ CoefficientsSet1,
uint32_t FilterCoefSet2RegBase,
CFilterCoefficientsNet^ CoefficientsSet2,
uint32_t FilterInfoRegBase,
CFilterPropertyNet^ FilterProp )
```

11.28.2.6 SetFilterParameterPermanent() [2/2] void SetFilterParameterPermanent (

```
uint32_t FilterCoefSet1RegBase,
uint32_t FilterCoefSet1DmaReg,
uint32_t FilterCoefSet2RegBase,
uint32_t FilterCoefSet2DmaReg,
uint32_t FilterInfoRegBase,
uint32_t FilterInfoDmaReg,
uint32_t EEPROMBase,
uint32_t EEPROMSize)
```

11.29 CFilterPropertyNet Class Reference

Public Member Functions

- CFilterPropertyNet (uint32_t CornerFrequenzymHz, uint32_t Order, FilterBandEnumNet FilterBand, FilterFamilyEnumNet FilterFamily, FilterTypeEnumNet FilterType, bool Active)
- ∼CFilterPropertyNet ()
- virtual System::String ^ ToString () override

Properties

- uint32_t CornerFrequencymHz [get]
- double CornerFrequency [get]
- uint32_t Order [get]
- FilterBandEnumNet FilterBand [get]
- FilterFamilyEnumNet FilterFamily [get]
- FilterTypeEnumNet FilterType [get]
- bool FilterActive [get]

11.29.1 Constructor & Destructor Documentation

11.29.1.1 CFilterPropertyNet() CFilterPropertyNet (


```
uint32_t CornerFrequenzymHz,
uint32_t Order,
FilterBandEnumNet FilterBand,
FilterFamilyEnumNet FilterFamily,
FilterTypeEnumNet FilterType,
bool Active )
```

```
11.29.1.2 ~CFilterPropertyNet() ~CFilterPropertyNet ( )
11.29.2 Member Function Documentation
11.29.2.1 ToString() virtual System::String ^ ToString ( ) [override], [virtual]
11.29.3 Property Documentation
11.29.3.1 CornerFrequency double CornerFrequency [get]
11.29.3.2 CornerFrequencymHz uint32_t CornerFrequencymHz [get]
11.29.3.3 FilterActive bool FilterActive [get]
11.29.3.4 FilterBand FilterBandEnumNet FilterBand [get]
11.29.3.5 FilterFamily FilterFamilyEnumNet FilterFamily [get]
11.29.3.6 FilterType FilterTypeEnumNet FilterType [get]
11.29.3.7 Order uint32_t Order [get]
```

11.30 CFluidControlDeviceNet Class Reference

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

Inheritance diagram for CFluidControlDeviceNet:

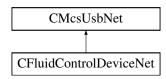

Public Member Functions

CFluidControlDeviceNet ()

Initialize a new instance of the CFluidControlDeviceNet class.

∼CFluidControlDeviceNet ()

Default destructor.

void SetValve (unsigned int value)

Open or Close valves.

void SetSingleValve (unsigned short valve, unsigned short onoff)

Opens or Closes a valve.

void SetDigout (unsigned int value)

Define the pattern on the Digital Output.

void SetPWM (unsigned int channel, unsigned int value)

Sets the duty cycle of the PWM output.

· void CalibrateThermocouple (unsigned int channel)

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)

Sets the proportinal constant for the Thermocouple.

• unsigned int GetValve ()

Gets the state of the valves.

unsigned short GetSingleValve (unsigned short valve)

Gets the state of a valve.

• unsigned int GetDigout ()

Gets the state of the digital output.

unsigned int GetPWM (unsigned int channel)

Gets the state of the PWM output.

• unsigned int GetAdc (unsigned int channel)

Reads an ADC Value.

• unsigned int GetDigin ()

Reads the digital input.

• int GetThermocoupleTemperature (unsigned int channel)

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

• int GetReferenceTemperature (unsigned int channel)

Reads the reference temperature for the Thermocouple.

unsigned int GetThermocoupleCalibration (unsigned int channel)

Gets the calibration constant for the Thermocouple ADC.

unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)

Reads the proportional constant for the Thermocouple.

Properties

• CMcsBus_VoltageModeNet^ McsBus_VoltageMode [get]

Additional Inherited Members

11.30.1 Detailed Description

CFluidControlDeviceNet is the class to control MCS FluidControl (FCB and FCX) device.

11.30.2 Constructor & Destructor Documentation

```
11.30.2.1 CFluidControlDeviceNet() CFluidControlDeviceNet ()
```

Initialize a new instance of the CFluidControlDeviceNet class.

```
11.30.2.2 ~CFluidControlDeviceNet() ~CFluidControlDeviceNet ()
```

Default destructor.

11.30.3 Member Function Documentation

```
11.30.3.1 CalibrateThermocouple() void CalibrateThermocouple ( unsigned int channel )
```

Calibrates the ADC which is used for the Thermocouple. For the calibration, Short circuit the Thermocouple and use this function to correct a possible offset of the ADC which measures the thermocouple.

Parameters

```
channel Thermocouple channel number.
```

```
11.30.3.2 GetAdc() unsigned int GetAdc ( unsigned int channel)
```

Reads an ADC Value.

channel	The ADC channel number to query.
---------	----------------------------------

Returns

The current ADC value.

```
11.30.3.3 GetDigin() unsigned int GetDigin ()
```

Reads the digital input.

Returns

The bit pattern of the state of the digital inputs.

```
11.30.3.4 GetDigout() unsigned int GetDigout ( )
```

Gets the state of the digital output.

Returns

The current state of the digital outputs as a bit pattern.

```
11.30.3.5 GetPWM() unsigned int GetPWM ( unsigned int channel)
```

Gets the state of the PWM output.

Returns

The current state of the PWM outputs duty cycle in permille.

```
11.30.3.6 GetReferenceTemperature() int GetReferenceTemperature ( unsigned int channel )
```

Reads the reference temperature for the Thermocouple.

channel	Thermocouple channel number.
---------	------------------------------

Returns

The temperature from the Thermocouple in 1/100 °C.

11.30.3.7 **GetSingleValve()** unsigned short GetSingleValve (unsigned short *valve*)

Gets the state of a valve.

Parameters

valve	number of valve
· a · · ·	mannoon on varvo

Returns

state of the valve

11.30.3.8 **GetThermocoupleCalibration()** unsigned int GetThermocoupleCalibration (unsigned int *channel*)

Gets the calibration constant for the Thermocouple ADC.

Parameters

channel	Thermocouple channel number.

Returns

The calibration constant for the Thermocouple ADC.

11.30.3.9 **GetThermocoupleNanovoltPerKelvin()** unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int *channel*)

Reads the proportional constant for the Thermocouple.

Returns

The proportional constant in Nanovolt per Kelvin.

11.30.3.10 **GetThermocoupleTemperature()** int GetThermocoupleTemperature (unsigned int *channel*)

Reads the temperature from Thermocouple. The functions gives the temperature difference between both Thermocouple junctions. To get the absolute temperature, add the reference temperature.

Parameters

channel	Thermocouple channel number.
---------	------------------------------

Returns

The temperature difference between both Thermocouple junctions in 1/100 ℃.

11.30.3.11 GetValve() unsigned int GetValve ()

Gets the state of the valves.

Returns

The current state of the valves as a bit pattern.

```
11.30.3.12 SetDigout() void SetDigout (
unsigned int value)
```

Define the pattern on the Digital Output.

Parameters

```
value bit pattern on the digital output.
```

```
11.30.3.13 SetPWM() void SetPWM (
unsigned int channel,
unsigned int value)
```

Sets the duty cycle of the PWM output.

channel	PWM channel number.
value	duty cycle of the PWM output in permille.

```
11.30.3.14 SetSingleValve() void SetSingleValve (
    unsigned short valve,
    unsigned short onoff)
```

Opens or Closes a valve.

Parameters

valve	number of valve to be changed.
-------	--------------------------------

Parameters

onoff open or close the valve.

```
11.30.3.15 SetThermocoupleNanovoltPerKelvin() void SetThermocoupleNanovoltPerKelvin (
    unsigned int channel,
    unsigned int value)
```

Sets the proportinal constant for the Thermocouple.

Parameters

channel	Thermocouple channel number.	
value	proportinal constant for the Thermocouple in Nanovolt per Kelvin.	

```
11.30.3.16 SetValve() void SetValve (
unsigned int value)
```

Open or Close valves.

value	bit pattern of valves which should be open.


11.30.4 Property Documentation

11.30.4.1 McsBus_VoltageMode CMcsBus_VoltageModeNet^ McsBus_VoltageMode [get]

11.31 CFYIDeviceNet Class Reference

CFYIDeviceNet is the class to control the MCS FYI device

Inheritance diagram for CFYIDeviceNet:

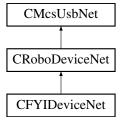

Public Member Functions

CFYIDeviceNet (void)

Properties

- CRobo_FYITemp_FunctionNet^ FYITemp [get]
- CRobo_FYIProgram_FunctionNet^ FYIProgram [get]
- CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.31.1 Detailed Description

CFYIDeviceNet is the class to control the MCS FYI device

11.31.2 Constructor & Destructor Documentation

```
11.31.2.1 CFYIDeviceNet() CFYIDeviceNet (
```

11.31.3 Property Documentation

11.31.3.1 FYIProgram CRobo_FYIProgram_FunctionNet^ FYIProgram [get]


11.31.3.2 FYITemp CRobo_FYITemp_FunctionNet^ FYITemp [get]

11.31.3.3 Sensor CMcsBus_SensorNet^ Sensor [get]

11.32 CGenericDevelopDeviceNet Class Reference

CGenericDevelopDeviceNet is the class to use during development of a new device.

Inheritance diagram for CGenericDevelopDeviceNet:

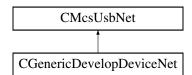

Public Member Functions

- CGenericDevelopDeviceNet (void)
 - Initialize a new instance of the CGenericDevelopDeviceNet class.
- ~CGenericDevelopDeviceNet (void)
- void SetValue (uint16_t value, uint16_t index)

Sets .

Parameters

value The value of the request.

Parameters

index The index of the request.

template < typename C > void SetBuffer (uint16_t value, uint16_t index, array < C >^ buffer)

void SetUByteBuffer (uint16_t value, uint16_t index, array< unsigned char >^ buffer)
 Sends an array of type unsigned char to the device.

Parameters

value The value of the reques	t.
-------------------------------	----

Parameters

index The index of the	he request.
------------------------	-------------

Parameters

buffer	The buffer to send.

• void SetByteBuffer (uint16_t value, uint16_t Index, array< char >^ buffer)

Sends an array of type char to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

Index	The index of the request.
mack	The mack of the reguest.

Parameters

ffer The buffer to :	send.
----------------------	-------

• void SetUShortBuffer (uint16_t value, uint16_t index, array< unsigned short >^ buffer)

Sends an array of type unsigned short to the device.

Parameters

value	The value of the request.
-------	---------------------------

index	The index of the request.
-------	---------------------------

• void SetShortBuffer (uint16_t value, uint16_t index, array< short $>^{\land}$ buffer) Sends an array of type short to the device.

Parameters

,	value	The value of the request.
---	-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

iffer The buffer to send	d.
--------------------------	----

• void SetUIntBuffer (uint16_t value, uint16_t index, array< unsigned int >^ buffer)

Sends an array of unsigned int to the device.

Parameters

value	The value of the request.

Parameters

index	The index of the request.
-------	---------------------------

Parameters

• void SetIntBuffer (uint16_t value, uint16_t index, array< int $>^{\land}$ buffer) Sends an array of type int to the device.

Parameters

index	The index of the request.
-------	---------------------------

Parameters

- template<typename C >
 array< C > ^ GetBuffer (uint16_t value, uint16_t index, int size)
- array< unsigned char > ^ GetUByteBuffer (uint16_t value, uint16_t index, int size)
 Gets an array of type unsigned char from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index The index of	of the request.
--------------------	-----------------


Parameters

size The size of the array

Returns

The array of data from the device.

array< char > ^ GetByteBuffer (uint16_t value, uint16_t index, int size)
 Gets an array of type char from the device.

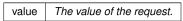

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

array< unsigned short > ^ GetUShortBuffer (uint16_t value, uint16_t index, int size)
 Gets an array of type unsigned short from the device.

Parameters

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

array< short > ^ GetShortBuffer (uint16_t value, uint16_t index, int size)
 Gets an array of type short from the device.

Parameters

value	The value of the request.
-------	---------------------------

index	The index of the request.
-------	---------------------------

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

- array< unsigned int > ^ GetUIntBuffer (uint16_t value, uint16_t index, int size)

Gets an array of type unsigned int from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index 7	The index of the request.
---------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

• array < int > $^{\land}$ GetIntBuffer (uint16_t value, uint16_t index, int size) Gets an array of type int from the device.

Parameters

value The value of the	request.
------------------------	----------

Parameters

index	The index of the request.
-------	---------------------------

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

• template<typename C>

void VendorOutRequest (uint8_t request, uint16_t value, uint16_t index, array< C >^ buffer)

• template<typename C >

array< C > ^ VendorInRequest (uint8 t request, uint16 t value, uint16 t index, int size)

- array< uint8_t > $^{\land}$ FindEndpoints (uint8_t type, uint8_t direction)
- IntPtr OpenPipe (uint8_t endpointAddress)

Open a Pipe to an USB Endpoint.

Parameters

endpointAddress Th	e Endpoint Number.
--------------------	--------------------

Returns

A handle to the endpoint.

• void ClosePipe (IntPtr pipeHandle)

Close a Pipe to an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.
------------	----------------------

• void ResetPipe (IntPtr pipeHandle)

Reset a Pipe to an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.
DIDENATION	i i ile elluboli il lialiule.

 template < typename C > array < C > ^ ReadPipe (IntPtr pipeHandle, uint32_t size)

Read data from an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.

Parameters

size Number of eler	ments to read.
---------------------	----------------

Returns

An array of data read.

template<typename C >
 void WritePipe (IntPtr pipeHandle, array< C >^ buffer)

Write data to an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.
------------	----------------------

Parameters

buffer An an	ray of data to write.
--------------	-----------------------

Additional Inherited Members

11.32.1 Detailed Description

CGenericDevelopDeviceNet is the class to use during development of a new device.

11.32.2 Constructor & Destructor Documentation

```
11.32.2.1 CGenericDevelopDeviceNet() CGenericDevelopDeviceNet (
void )
```

Initialize a new instance of the CGenericDevelopDeviceNet class.

```
11.32.2.2 ~CGenericDevelopDeviceNet() ~CGenericDevelopDeviceNet (
```

11.32.3 Member Function Documentation

Close a Pipe to an USB Endpoint.

pipeHandle	The endpoint handle.

Gets an array of type char from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

	index	The	index	of the	request.
--	-------	-----	-------	--------	----------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

Gets an array of type int from the device.

Parameters

value The value	of the request.
-----------------	-----------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

Gets an array of type short from the device.

Parameters

value	The value of the request.

index The index of the reques

Returns

The array of data from the device.

Gets an array of type unsigned char from the device.

Parameters

value The value of th	e request.
-----------------------	------------

Parameters

index The index of the req	uest.
----------------------------	-------

Parameters

Returns

The array of data from the device.

Gets an array of type unsigned int from the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

index	The index of the request.
-------	---------------------------

Parameters

size	The size of the array.
------	------------------------

Returns

The array of data from the device.

Gets an array of type unsigned short from the device.

Parameters

1	The second secon
vaiue	The value of the request.

index The index of the reque

size The size of the array.	' .
-----------------------------	------------

Returns

The array of data from the device.

Open a Pipe to an USB Endpoint.

Parameters

endpointAddress	The Endpoint Number.
-----------------	----------------------

Returns

A handle to the endpoint.

Read data from an USB Endpoint.

Parameters

pipeHandle	The endpoint handle.

Parameters

size Number of elements to read.

Returns

An array of data read.

Reset a Pipe to an USB Endpoint.

Parameters

Sends an array of type char to the device.

Parameters

value The value of the request.

Index The index of	the request.
--------------------	--------------

buffer The buffer to send.

Sends an array of type int to the device.

Parameters

value	The value of the request.
-------	---------------------------

Parameters

Parameters

buffer The buffer to send.

Sends an array of type short to the device.

value	The value of the request.
-------	---------------------------

ir	idex	The index of the request.
----	------	---------------------------

Parameters

```
buffer The buffer to send.
```

Sends an array of type unsigned char to the device.

Parameters

Parameters

Parameters

11.32.3.18 SetUIntBuffer() void SetUIntBuffer (uint16_t value,

```
uint16_t index,
array< unsigned int >^{\land} buffer)
```

Sends an array of unsigned int to the device.

Parameters

value The value of the request.	
---------------------------------	--

Parameters

Parameters

buffer	The buffer to send.
--------	---------------------

Sends an array of type unsigned short to the device.

Parameters

value	The value of the request.
-------	---------------------------

dex The index of the request.	index
-------------------------------	-------

buffer The buffer to send.

Sets.

Parameters

The value of the request.	value
---------------------------	-------

Parameters

```
index The index of the request.
```

```
uint16_t value,
uint16_t index,
array< C >^ buffer )
```

```
11.32.3.23 WritePipe() void WritePipe (

IntPtr pipeHandle,

array < C >^ buffer)
```

Write data to an USB Endpoint.

pipeHandle	The endpoint handle.
------------	----------------------


Parameters

buffer	An array of data to write.
--------	----------------------------

11.33 CGilsonDeviceNet Class Reference

CGilsonDeviceNet is the class to control a Gilson device.

Inheritance diagram for CGilsonDeviceNet:

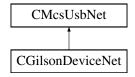

Public Member Functions

- CGilsonDeviceNet (void)
 - Initialize a new instance of the CGilsonDeviceNet class.
- ∼CGilsonDeviceNet (void)
- void ConnectSlave (byte ID)
- void SendImmediate (wchar_t command)
- String \(^\) SendImmediateGetResponse (wchar_t command)
- void SendBuffered (String[^] command)
- String ^ GetLastAnswer ()

Protected Attributes

• CGilsonDevice * m_pGilsonDevice

Additional Inherited Members

11.33.1 Detailed Description

CGilsonDeviceNet is the class to control a Gilson device.

11.33.2 Constructor & Destructor Documentation

```
11.33.2.1 CGilsonDeviceNet() CGilsonDeviceNet (
void )
```

Initialize a new instance of the CGilsonDeviceNet class.

```
11.33.2.2 ~CGilsonDeviceNet() ~CGilsonDeviceNet ( void )
```

11.33.3 Member Function Documentation

```
11.33.3.1 ConnectSlave() void ConnectSlave (
byte ID )
```

```
11.33.3.2 GetLastAnswer() String ^ GetLastAnswer ( )
```

```
11.33.3.3 SendBuffered() void SendBuffered ( String^{\wedge} command )
```


```
11.33.3.5 SendImmediateGetResponse() String ^ SendImmediateGetResponse ( wchar_t command )
```

11.33.4 Member Data Documentation

11.33.4.1 m_pGilsonDevice CGilsonDevice* m_pGilsonDevice [protected]

11.34 CGrapheneASICDeviceNet Class Reference

Inheritance diagram for CGrapheneASICDeviceNet:

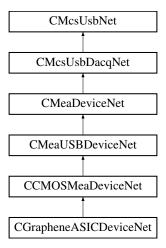

Public Member Functions

- CGrapheneASICDeviceNet (void)
- ∼CGrapheneASICDeviceNet ()
- void SetBaseSamplerate (int BaseSamplerate) override
- array< int > $^{\wedge}$ GetAvailableBaseSamplerates () override
- void SetRegionOfInterests (System::Collections::Generic::Dictionary < int, CCMOSMeaDeviceNet::CRegionOfInterestRect[^] > rois) override

Additional Inherited Members

11.34.1 Constructor & Destructor Documentation

```
11.34.1.1 CGrapheneASICDeviceNet() CGrapheneASICDeviceNet (
```

11.34.1.2 ~CGrapheneASICDeviceNet() ~CGrapheneASICDeviceNet ()

11.34.2 Member Function Documentation

```
11.34.2.1 GetAvailableBaseSamplerates() array<int> ^ GetAvailableBaseSamplerates ( ) [override], [virtual]
```


Reimplemented from CCMOSMeaDeviceNet.

Reimplemented from CCMOSMeaDeviceNet.

11.35 CGrapheneFunctionNet Class Reference

CGrapheneFunctionNet is the class to control Graphene device functions

Inheritance diagram for CGrapheneFunctionNet:

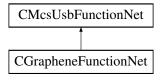

Public Member Functions

• CGrapheneFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] pGraphene ← FunctionPointerContainer)

Initializes a new instance of the CGrapheneFunctionNet class.

- CGrapheneFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CGrapheneFunctionNet ()
- !CGrapheneFunctionNet ()
- void GetVdVsDAC ([System::Runtime::InteropServices::Out]int16_t% Vd, [System::Runtime::Interop← Services::Out]int16_t% Vs)

Gets Vd and Vs

void GetVdVsDAC (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% Vd, [System::←
Runtime::InteropServices::Out]int16_t% Vs)

Gets Vd and Vs

void SetVdVsDAC (int16 t Vd, int16 t Vs)

Sets Vd and Vs

void SetVdVsDAC (uint32 t Headstage, int16 t Vd, int16 t Vs)

Sets Vd and VS

• bool GetVoltageReached ()

Gets the reached voltage

 bool GetVoltageReached (uint32_t Headstage) Gets the reached voltage • int32_t GetVoltageRange () Gets the voltage range int32_t GetVoltageRange (uint32_t Headstage) Gets the voltage range void SetVoltageRange (int32 t range) Sets the voltage range void SetVoltageRange (uint32_t Headstage, int32_t range) Sets the voltage range int32_t GetVoltageResolution () Gets the voltage resolution int32_t GetVoltageResolution (uint32_t Headstage) Gets the voltage resolution void SetVoltageResolution (int32 t resolution) Sets the voltage resolution void SetVoltageResolution (uint32_t Headstage, int32_t resolution) Sets the voltage resolution void GetDACOffset ([System::Runtime::InteropServices::Out]int16 t% offset vd, [System::Runtime::
 — InteropServices::Out]int16_t% offset_vs) Gets the DAC offset • void GetDACOffset (uint32_t Headstage, [System::Runtime::InteropServices::Out]int16_t% offset_vd, [System::Runtime::InteropServices::Out]int16_t% offset_vs) Gets the DAC offset void SetDACOffset (int16_t offset_vd, int16_t offset_vs) Sets the DAC offset void SetDACOffset (uint32_t Headstage, int16_t offset_vd, int16_t offset_vs) Set the DAC offset void GetVdVs ([System::Runtime::InteropServices::Out]int32_t% Vd, [System::Runtime::InteropServices:: Outlint32 t% Vs) Gets Vd and Vs void GetVdVs (uint32 t Headstage, [System::Runtime::InteropServices::Out]int32 t% Vd, [System::← Runtime::InteropServices::Out]int32 t% Vs) Gets Vd and Vs void SetVdVs (int32 t Vd, int32 t Vs) Sets Vd and Vs void SetVdVs (uint32_t Headstage, int32_t Vd, int32_t Vs) Sets Vd and Vs void SetVdsVgs (int32_t Vds, int32_t Vgs) Sets Vds and Vgs void SetVdsVgs (uint32_t Headstage, int32_t Vds, int32_t Vgs) Sets Vds and Vqs int32 t GetCur2VolResistance () Gets the resistance of the current to voltage converter int32_t GetCur2VolResistance (uint32_t Headstage)

Gets the resistance of the current to voltage converter

• int32_t GetADCOffset (uint32_t Channel)

Gets the ADC offset

int32_t GetADCOffset (uint32_t Headstage, uint32_t Channel)

Gets the ADC offset

void SetADCOffset (uint32_t Channel, int32_t offset)

```
Sets the ADC offset
```

void SetADCOffset (uint32_t Headstage, uint32_t Channel, int32_t offset)

Sets the ADC offset

• int32_t GetADCGain (uint32_t Channel)

Gets the ADC offset

int32_t GetADCGain (uint32_t Headstage, uint32_t Channel)

Gets the ADC offset

void SetADCGain (uint32_t Channel, int32_t gain)

Sets the ADC offset

· void SetADCGain (uint32 t Headstage, uint32 t Channel, int32 t gain)

Sets the ADC offset

void GetVdsVgs ([System::Runtime::InteropServices::Out]int32_t% Vds, [System::Runtime::Interop←
 Services::Out]int32_t% Vgs)

Gets Vds and Vgs

 void GetVdsVgs (uint32_t Headstage, [System::Runtime::InteropServices::Out]int32_t% Vds, [System::← Runtime::InteropServices::Out]int32_t% Vqs)

Gets Vds and Vgs

Additional Inherited Members

11.35.1 Detailed Description

CGrapheneFunctionNet is the class to control Graphene device functions

11.35.2 Constructor & Destructor Documentation

```
11.35.2.1 CGrapheneFunctionNet() [1/2] CGrapheneFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pGrapheneFunctionPointerContainer)
```

Initializes a new instance of the CGrapheneFunctionNet class.

```
11.35.2.2 CGrapheneFunctionNet() [2/2] CGrapheneFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.35.2.3 ~CGrapheneFunctionNet() virtual ~CGrapheneFunctionNet ( ) [virtual]
```

```
11.35.2.4 "!CGrapheneFunctionNet() !CGrapheneFunctionNet ( )
```

11.35.3 Member Function Documentation

```
11.35.3.1 GetADCGain() [1/2] int32_t GetADCGain ( uint32_t Channel )
```

Gets the ADC offset

Channel	The channel to query.
---------	-----------------------

Returns

The gain

Gets the ADC offset

Parameters

Headstage	The headstage to query.
Channel	The channel to query.

Returns

The gain

11.35.3.3 GetADCOffset() [1/2] int32_t GetADCOffset (uint32_t Channel)

Gets the ADC offset

Parameters

Channel	The channel to query.
---------	-----------------------

Returns

The offset

Gets the ADC offset

Headstage	The headstage to query.
Channel	The channel to query.

Returns

The offset

11.35.3.5 GetCur2VolResistance() [1/2] int32_t GetCur2VolResistance ()

Gets the resistance of the current to voltage converter

Returns

The resistance in Ohm

11.35.3.6 GetCur2VolResistance() [2/2] int32_t GetCur2VolResistance (uint32_t Headstage)

Gets the resistance of the current to voltage converter

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The resistance in Ohm

Gets the DAC offset

offset_vd	Vd offset in DAC Units
offset_vs	Vs offset in DAC Units

11.35.3.8 GetDACOffset() [2/2] void GetDACOffset (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int16_t% offset_vd,
[System::Runtime::InteropServices::Out] int16_t% offset_vs )
```

Gets the DAC offset

Parameters

Headstage	The headstage to query.
offset_vd	Vd offset in DAC Units
offset_vs	Vs offset in DAC Units

11.35.3.9 GetVdsVgs() [1/2] void GetVdsVgs (

```
[System::Runtime::InteropServices::Out] int32_t% Vds, [System::Runtime::InteropServices::Out] int32_t% Vgs )
```

Gets Vds and Vgs

Parameters

Vds	Vds in μV
Vgs	Vgs in μV

11.35.3.10 GetVdsVgs() [2/2] void GetVdsVgs (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int32_t% Vds,
[System::Runtime::InteropServices::Out] int32_t% Vgs )
```

Gets Vds and Vgs

Parameters

Headstage	The headstage to query.
Vds	Vds in μV
Vgs	Vgs in μV

11.35.3.11 GetVdVs() [1/2] void GetVdVs (

```
[System::Runtime::InteropServices::Out] int32_t% Vd, [System::Runtime::InteropServices::Out] int32_t% Vs)
```

Gets Vd and Vs

Vd	Vd in μV
Vs	Vs in μV

11.35.3.12 GetVdVs() [2/2] void GetVdVs (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int32_t% Vd,
[System::Runtime::InteropServices::Out] int32_t% Vs )
```

Gets Vd and Vs

Parameters

Headstage	The headstage to query.
Vd	Vd in μV
Vs	Vs in μV

11.35.3.13 GetVdVsDAC() [1/2] void GetVdVsDAC (

```
[System::Runtime::InteropServices::Out] int16_t% Vd, [System::Runtime::InteropServices::Out] int16_t% Vs)
```

Gets Vd and Vs

Parameters

Vd	Vd in DAC Units
Vs	Vs in DAC Units

11.35.3.14 GetVdVsDAC() [2/2] void GetVdVsDAC (

```
uint32_t Headstage,
[System::Runtime::InteropServices::Out] int16_t% Vd,
[System::Runtime::InteropServices::Out] int16_t% Vs )
```

Gets Vd and Vs

Headstage	The headstage to query.
Vd	Vd in DAC Units
Vs	Vs in DAC Units

```
11.35.3.15 GetVoltageRange() [1/2] int32_t GetVoltageRange ( )
```

Gets the voltage range

Returns

The voltage range in mV

Gets the voltage range

Parameters

eadstage The headstage to query.	Headstage
----------------------------------	-----------

Returns

The voltage range in mV

11.35.3.17 GetVoltageReached() [1/2] bool GetVoltageReached ()

Gets the reached voltage

Returns

the reached voltage

```
11.35.3.18 GetVoltageReached() [2/2] bool GetVoltageReached ( uint32_t Headstage )
```

Gets the reached voltage

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The reached voltage

11.35.3.19 GetVoltageResolution() [1/2] int32_t GetVoltageResolution ()

Gets the voltage resolution

Returns

The voltage resolution in $\mu V/\text{digit}$

11.35.3.20 GetVoltageResolution() [2/2] int32_t GetVoltageResolution (uint32_t Headstage)

Gets the voltage resolution

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The voltage resolution in $\mu V/\text{digit}$

11.35.3.21 SetADCGain() [1/2] void SetADCGain (uint32_t Channel, int32_t gain)

Sets the ADC offset

Parameters

Channel	The channel to query.
gain	The gain

Sets the ADC offset

Headstage	The headstage to query.
Channel	The channel to query.
gain	The gain

Sets the ADC offset

Parameters

Channel	The channel to query.
offset	The offset

Sets the ADC offset

Parameters

Headstage	The headstage to query.
Channel	The channel to query.
offset	The offset

Sets the DAC offset

Parameters

offset_vd	Vd offset in DAC Units
offset_vs	Vs offset in DAC Units

Set the DAC offset

Headstage	The headstage to query.
offset_vd	Vd offset in DAC Units
offset_vs	Vs offset in DAC Units

Sets Vds and Vgs

Parameters

Vds	Vds in μV
Vgs	Vgs in μV

Sets Vds and Vgs

Parameters

Headstage	The headstage to query.
Vds	Vds in μV
Vgs	Vgs in μV

Sets Vd and Vs

Vd	Vd in μV
Vs	Vs in μV

int32_t *Vs*)

Sets Vd and Vs

Parameters

Headstage	The headstage to query.
Vd	Vd in μV
Vs	Vs in μV

Sets Vd and Vs

Parameters

Vd	Vd in DAC Units
Vs	Vs in DAC Units

Sets Vd and VS

Parameters

Headstage	The headstage to query.
Vd	Vd in DAC Units
Vs	Vs in DAC Units

```
11.35.3.33 SetVoltageRange() [1/2] void SetVoltageRange ( int32_t range )
```

Sets the voltage range

range	The voltage range in mV
-------	-------------------------

Sets the voltage range

Parameters

Headstage	The headstage to query.
range	The voltage range in mV

Sets the voltage resolution

Parameters

resolution	The voltage resolution in μV/digit
------------	------------------------------------

```
11.35.3.36 SetVoltageResolution() [2/2] void SetVoltageResolution ( uint32_t Headstage, int32_t resolution )
```

Sets the voltage resolution


Parameters

Headstage	The headstage to query.
resolution	The voltage resolution in μV/digit

11.36 CHiClampDeviceNet Class Reference

CHiClampDeviceNet is the to control the MCS HiClamp device

Inheritance diagram for CHiClampDeviceNet:

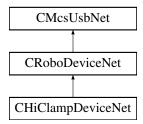

Public Member Functions

• CHiClampDeviceNet (void)

Properties

• CRoboDacqNet^ RoboDacq [get]

Additional Inherited Members

11.36.1 Detailed Description

CHiClampDeviceNet is the to control the MCS HiClamp device

11.36.2 Constructor & Destructor Documentation


```
11.36.2.1 CHiClampDeviceNet() CHiClampDeviceNet (
void )
```

11.36.3 Property Documentation

11.36.3.1 RoboDacq CRoboDacqNet^ RoboDacq [get]

11.37 CHLADacqNet Class Reference

Inheritance diagram for CHLADacqNet:

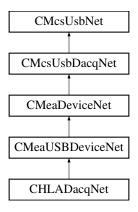

Public Member Functions

CHLADacqNet (void)

Additional Inherited Members


11.37.1 Constructor & Destructor Documentation

```
11.37.1.1 CHLADacqNet() CHLADacqNet ( void )
```

11.38 CHLADeviceNet Class Reference

CHLADeviceNet is the to control the MCS HLA device

Inheritance diagram for CHLADeviceNet:

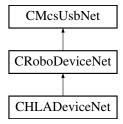

Public Member Functions

CHLADeviceNet (void)

Properties

- CHLADacqNet^ HLADacq [get]
- CSerialPortNet^ SerialPort [get]

Additional Inherited Members

11.38.1 Detailed Description

CHLADeviceNet is the to control the MCS HLA device

11.38.2 Constructor & Destructor Documentation

```
11.38.2.1 CHLADeviceNet() CHLADeviceNet (
```

11.38.3 Property Documentation

```
11.38.3.1 HLADacq CHLADacqNet^ HLADacq [get]
```

```
11.38.3.2 SerialPort CSerialPortNet^ SerialPort [get]
```

11.39 CMcsUsbDacqNet::CHWInfo Class Reference

Class to provide hardware information about the device.

Classes

· class CVoltageRangeInfoNet

Public Member Functions

- CHWInfo (CMcsUsbDacqNet[^] device)

Get the number of analog channels the device supports.

Get the number of digital channels the device supports.

virtual bool IsDigitalChannelDedicated ()

Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

- virtual uint32_t GetAvailableSampleRates ([System::Runtime::InteropServices::Out]System::Collections::
 Generic::List< int32_t >^% sampleRates)
- virtual System::Collections::Generic::List< int32_t > ^ GetAvailableVoltageRangesInMicroVolt (int milliGain)

 Gets a List of voltage ranges the device supports.
- virtual System::Collections::Generic::List< CVoltageRangeInfoNet[^]> [^] GetAvailableVoltageRangesInMicroVoltAndStringsInM (int milliGain)

Gets a List of voltage ranges the device supports.

11.39.1 Detailed Description

Class to provide hardware information about the device.

11.39.2 Constructor & Destructor Documentation

```
11.39.2.1 CHWInfo() CHWInfo (

CMcsUsbDacqNet^ device)
```

11.39.3 Member Function Documentation

Gets a List of voltage ranges the device supports.

The List is scaled by the gain factor given as argument to this function. Use "1000" as scale factor for backwards compatibility. To get a list of voltage ranges for the headstage, obtain the gain of the headstage with the Get← Gain() call and use the result in the milliGain parameter. To get a list of voltage ranges for the analog inputs of the interfaceboard, obtain the gain of the analog inputs with the GetAnalogGain() call and use the result in the milliGain parameter.

Parameters

milliGain The gain factor (in milliGain) used to scale the list of voltage ranges.

Returns

List of voltage ranges in µV.

Gets a List of voltage ranges the device supports.

The List is scaled by the gain factor given as argument to this function. Use "1000" as scale factor for backwards compatibility. Each list entry contains the voltage range as an integer and as a string. To get a list of voltage ranges for the headstage, obtain the gain of the headstage with the GetGain() call and use the result in the milliGain parameter. To get a list of voltage ranges for the analog inputs of the interfaceboard, obtain the gain of the analog inputs with the GetAnalogGain() call and use the result in the milliGain parameter.

milliGain The gain factor (in milliGain) used to scale the list of voltage ranges.
--

Returns

List of voltage ranges in μV .

Get the number of analog channels the device supports.

Parameters

numberOfChannels	Number of analog channels the device supports.
------------------	--

Returns

Error Status. 0 on success.

11.39.3.5 **GetNumberOfHWDigitalChannels()** virtual uint32_t GetNumberOfHWDigitalChannels ([System::Runtime::InteropServices::Out] int% numberOfChannels) [virtual]

Get the number of digital channels the device supports.

Parameters

numberOfChannels	Number of digital channels the device supports.
number Of Charmers	Number of digital charmers the device supports.

Returns

Error Status. 0 on success.

11.39.3.6 IsDigitalChannelDedicated() virtual bool IsDigitalChannelDedicated () [virtual]


Query if the digital channel replaces an analog channel when enabled (e.g. on MC_Card) or adds a channel link on USB devices.

Returns

false when the digital channel replaces an analog channel when enabled, true when the digital channels is appended to the analog channels when enabled.

11.40 CIntanMea FunctionNet Class Reference

Inheritance diagram for CIntanMea_FunctionNet:

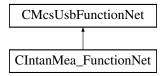

Public Member Functions

- CIntanMea_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] intalMea_Function
 — PointerContainer)
- CIntanMea_FunctionNet (CMcsUsbNet[^] mcsusb)
- int GetUpperFrequencyByIndex (unsigned short index)
- int GetLowerFrequencyByIndex (unsigned short index)
- int64_t GetDSPHighPassByIndex (unsigned short index)
- int GetIntanRegister (unsigned short chip, unsigned short registernumber)
- int GetImpedanceResult (unsigned short channel)
- void SetBandwidthByIndex (int upper_index, int lower_index)
- void SetDSPHighPassByIndex (int index)
- void AmplifierSettle ()
- void SetIntanRegister (unsigned short register_number, int value)
- void SetDiagnosticMode (unsigned char onoff)
- void BeginImpedanceCheck (array< int >^ config_values)

Additional Inherited Members


11.40.1 Constructor & Destructor Documentation

```
11.40.1.2 CIntanMea_FunctionNet() [2/2] CIntanMea_FunctionNet (
             CMcsUsbNet^ mcsusb )
11.40.2 Member Function Documentation
11.40.2.1 AmplifierSettle() void AmplifierSettle ( )
11.40.2.2 BeginImpedanceCheck() void BeginImpedanceCheck (
             array < int >^{\land} config\_values)
11.40.2.3 GetDSPHighPassByIndex() int64_t GetDSPHighPassByIndex (
             unsigned short index )
11.40.2.4 GetImpedanceResult() int GetImpedanceResult (
             unsigned short channel )
11.40.2.5 GetIntanRegister() int GetIntanRegister (
             unsigned short chip,
             unsigned short registernumber )
11.40.2.6 GetLowerFrequencyByIndex() int GetLowerFrequencyByIndex (
             unsigned short index)
11.40.2.7 GetUpperFrequencyByIndex() int GetUpperFrequencyByIndex (
             unsigned short index)
11.40.2.8 SetBandwidthByIndex() void SetBandwidthByIndex (
             int upper_index,
             int lower_index )
```

11.41 CInterfaceboard2FunctionNet Class Reference

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboard2FunctionNet:

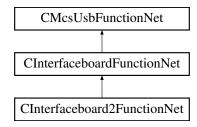

Public Member Functions

 CInterfaceboard2FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pInterfaceboard2← FunctionPointerContainer)

Initializes a new instance of the CInterfaceboard2FunctionNet class.

- CInterfaceboard2FunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CInterfaceboard2FunctionNet ()
- !CInterfaceboard2FunctionNet ()
- void SetloVoltage (IoVoltageEnumNet ioVoltage)

Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

loVoltageEnumNet GetloVoltage ()

Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.

Additional Inherited Members

11.41.1 Detailed Description

CInterfaceboard2FunctionNet is the class to control the Interfaceboard

Parameters ioVoltage

Enum for the I/O Voltage (3.3V or 5.0V).


11.41.2 Constructor & Destructor Documentation

```
11.41.2.1 CInterfaceboard2FunctionNet() [1/2] CInterfaceboard2FunctionNet (
              CMcsUsbNet^ mcsusb,
              {\tt CMcsUsbFunctionPointerContainer}^{\wedge} \ pInterfaceboard2FunctionPointerContainer \ )
Initializes a new instance of the CInterfaceboard2FunctionNet class.
11.41.2.2 CInterfaceboard2FunctionNet() [2/2] CInterfaceboard2FunctionNet (
              CMcsUsbNet^ mcsusb )
11.41.2.3 ~CInterfaceboard2FunctionNet() virtual ~CInterfaceboard2FunctionNet ( ) [virtual]
11.41.2.4 "!CInterfaceboard2FunctionNet() !CInterfaceboard2FunctionNet ( )
11.41.3 Member Function Documentation
11.41.3.1 GetloVoltage() IoVoltageEnumNet GetIoVoltage ( )
Gets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.
Returns
     Enum for the IO Voltage (3.3V or 5.0V).
11.41.3.2 SetloVoltage() void SetIoVoltage (
              IoVoltageEnumNet ioVoltage )
Sets the I/O Voltage level for the IFB2 digital and AUX ports, default is 3.3V.
```

11.42 CInterfaceboardFunctionNet Class Reference

CInterfaceboardFunctionNet is the class to control the Interfaceboard

Inheritance diagram for CInterfaceboardFunctionNet:

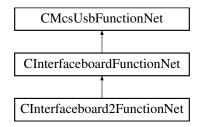

Public Member Functions

CInterfaceboardFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pInterfaceboard
 —
 FunctionPointerContainer)

Initializes a new instance of the CInterfaceboardFunctionNet class.

- CInterfaceboardFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CInterfaceboardFunctionNet ()
- !CInterfaceboardFunctionNet ()
- void SetCardinalDacqSamplerate (uint32_t samplerate)

Sets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz

uint32_t GetCardinalDacqSamplerate ()

Gets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz

void SetCardinalStgOutputrate (uint32_t outputrate)

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

• uint32_t GetCardinalStgOutputrate ()

Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz

Additional Inherited Members

11.42.1 Detailed Description

CInterfaceboardFunctionNet is the class to control the Interfaceboard

11.42.2 Constructor & Destructor Documentation

```
11.42.2.1 CInterfaceboardFunctionNet() [1/2] CInterfaceboardFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pInterfaceboardFunctionPointerContainer)
```

Initializes a new instance of the CInterfaceboardFunctionNet class.

```
11.42.2.2 CInterfaceboardFunctionNet() [2/2] CInterfaceboardFunctionNet (
              CMcsUsbNet^ mcsusb )
11.42.2.3 ~CInterfaceboardFunctionNet() virtual ~CInterfaceboardFunctionNet () [virtual]
11.42.2.4 "!CInterfaceboardFunctionNet() !CInterfaceboardFunctionNet ( )
11.42.3 Member Function Documentation
11.42.3.1 GetCardinalDacqSamplerate() uint32_t GetCardinalDacqSamplerate ( )
Gets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz
Returns
     The samplerate in Hz.
11.42.3.2 GetCardinalStgOutputrate() uint32_t GetCardinalStgOutputrate ( )
Gets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz
Returns
     The output rate in Hz.
11.42.3.3 SetCardinalDacqSamplerate() void SetCardinalDacqSamplerate (
              uint32_t samplerate )
Sets the fundamental/cardinal data aquisition samplerate of the Interfaceboard, default is 50 kHz
Parameters
              The samplerate in Hz.
 samplerate
```

11.42.3.4 SetCardinalStgOutputrate() void SetCardinalStgOutputrate (uint32_t outputrate)

Sets the fundamental/cardinal STG output rate of the Interfaceboard, default is 50 kHz


Parameters

outputrate The output rate in Hz.

11.43 CLIH3DeviceNet Class Reference

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

Inheritance diagram for CLIH3DeviceNet:

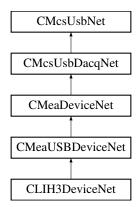

Public Member Functions

· CLIH3DeviceNet ()

Initializes a new instance of the CLIH3DeviceNet class.

- virtual ∼CLIH3DeviceNet ()
- !CLIH3DeviceNet ()
- void DummyCommand (uint32_t dummyParameter)

Dummy command to show how to use the DLL.

void SetEEpromPage (uint32_t EEpromStartAddress, array < int8_t >^ EEpromData, LIH30_EPC10_Bus_EnumNet epc10bus)

Writes into EEprom on the EPC10 EEPROM

array< int8_t > ^ GetEEpromPage (uint32_t EEpromStartAddress, int EEpromData_Length, LIH30_EPC10_Bus_EnumNet epc10bus)

Reads the requested amount of EEprom byte from the EPC10 EEPROM

void SetSampleInterval (uint32_t SampleInterval)

Sets the Sample Interval for the DACQ and Stimulation

• uint32 t GetSampleInterval ()

Gets the Sample Interval for the DACQ and Stimulation

void SetAdcOffset (LIH30_ADC_Channel_EnumNet AdcChannel, int32_t Offset)

Sets the ADC offset of the DACQ for a single channel

• int32_t GetAdcOffset (LIH30_ADC_Channel_EnumNet AdcChannel)

Gets the ADC offset of the DACQ for a single channel

• void SetAdcOffsetPermanent (LIH30_ADC_Channel_EnumNet AdcChannel)

Writes the ADC offset of the DACQ for a single channel to permanent EEProm memory

void ErasePermanentAdcOffset (LIH30 ADC Channel EnumNet AdcChannel)

Delets the ADC offset of the DACQ for a single channel in permanent EEProm memory

uint32_t ReadClipping (LIH30_EPC10_Bus_EnumNet epc10bus)

Gets the clipping information

void SetDigOutState (uint16_t DigOutState)

Writes to the LIH30 digital output

• uint16_t GetDigInState ()

Reads from the LIH30 digital input

• void SendCommand (LIH30_EPC10_Bus_EnumNet epc10bus, uint16_t Command)

Send command to the EPC10

uint16_t GetDacqRunStatus ()

Gets the data acquisition running status

void SetDacUseIdleValue (uint32_t DacChannel, bool UseIdle)

Sets if the DAC Idle value is used after stimulation

bool GetDacUseldleValue (uint32_t DacChannel)

Gets if the DAC Idle value is used after stimulation

void SetDacIdleValue (uint32 t DacChannel, int32 t IdleValue)

Sets the DAC Idle value

• int32 t GetDacIdleValue (uint32 t DacChannel)

Gets the DAC Idle value

void EnableUserTrigger (bool enable)

Enables the User Trigger

• bool IsUserTriggerEnabled ()

Is the User Trigger enabled

void SetDacOffset (LIH30_DAC_Channel_EnumNet DacChannel, int32_t Offset)

Sets the offset of a DAC channel.

int32_t GetDacOffset (LIH30_DAC_Channel_EnumNet DacChannel)

Gets the offset of a DAC channel.

void SetDacOffsetPermanent (LIH30_DAC_Channel_EnumNet DacChannel)

Writes the DAC offset of the STG for a single channel to permanent EEProm memory

void ErasePermanentDacOffset (LIH30_DAC_Channel_EnumNet DacChannel)

Delets the DAC offset of the STG for a single channel in permanent EEProm memory
• void SetAudioOutDacParameter (uint32_t Frequency, uint32_t Amplification)

Sets the parameter of the audio DAC output.

• void GetAudioOutDacParameter ([System::Runtime::InteropServices::Out]uint32_t% Frequency, [System::

Runtime::InteropServices::Out]uint32_t% Amplification)

Gets the parameter of the audio DAC output.

String ^ ReadUARTData ()

Reads the config string from the device.

void WriteUARTData (String[^] commandString)

Write the command string to the device.

Properties

CStimulusFunctionNet[^] StimulusFunction [get]

Additional Inherited Members

11.43.1 Detailed Description

CLIH3DeviceNet is the class to access the HEKA LIH3 device.

11.43.2 Constructor & Destructor Documentation

```
11.43.2.1 CLIH3DeviceNet() CLIH3DeviceNet ()
```

Initializes a new instance of the CLIH3DeviceNet class.

```
11.43.2.2 ~CLIH3DeviceNet() virtual ~CLIH3DeviceNet () [virtual]
```

```
11.43.2.3 "!CLIH3DeviceNet() !CLIH3DeviceNet ()
```

11.43.3 Member Function Documentation

```
11.43.3.1 DummyCommand() void DummyCommand ( uint32_t dummyParameter)
```

Dummy command to show how to use the DLL.

Parameters

dummyParameter | parameter to send to the device

```
11.43.3.2 EnableUserTrigger() void EnableUserTrigger ( bool enable)
```

Enables the User Trigger

Parameters

enable Enable

Delets the ADC offset of the DACQ for a single channel in permanent EEProm memory

Parameters

```
AdcChannel The ADC channel
```

Delets the DAC offset of the STG for a single channel in permanent EEProm memory

Parameters

DacChannel TI	he DAC channel
---------------	----------------

Gets the ADC offset of the DACQ for a single channel

Parameters

AdcChannel	The ADC channel

Returns

The offset for the given channel number

Gets the parameter of the audio DAC output.

Frequency	Frequency(1 - 25000 Hz)
Amplification	Amplification(0 - 0xFFFF)

```
11.43.3.7 GetDacIdleValue() int32_t GetDacIdleValue ( uint32_t DacChannel )
```

Gets the DAC Idle value

Parameters

```
DacChannel The DAC channel
```

Returns

The idle value

Gets the offset of a DAC channel.

Parameters

```
DacChannel The DAC channel
```

Returns

The offset for the given channel number

```
\textbf{11.43.3.9} \quad \textbf{GetDacqRunStatus()} \quad \texttt{uint16\_t GetDacqRunStatus ()}
```

Gets the data acquisition running status

Returns

The status (1: running / 0: stopped)

```
11.43.3.10 GetDacUseldleValue() bool GetDacUseIdleValue ( uint32_t DacChannel )
```

Gets if the DAC Idle value is used after stimulation

DacChannel The DAC channel

Returns

Use idle value

11.43.3.11 GetDigInState() uint16_t GetDigInState ()

Reads from the LIH30 digital input

Returns

The bit mask defining the digital input state

Reads the requested amount of EEprom byte from the EPC10 EEPROM

Parameters

EEpromStartAddress	start address of memory area to read from
EEpromData_Length	The maximal length of EEpromData.
epc10bus	The EPC10 bus

Returns

pointer to internal memory for the requested amount of data

11.43.3.13 GetSampleInterval() uint32_t GetSampleInterval ()

Gets the Sample Interval for the DACQ and Stimulation

Returns

Sample Interval configured on the device

11.43.3.14 IsUserTriggerEnabled() bool IsUserTriggerEnabled ()

Is the User Trigger enabled

Returns

Enabled

Gets the clipping information

Parameters

```
epc10bus The EPC10 bus
```

Returns

The clipping value

11.43.3.16 ReadUARTData() String $^{\wedge}$ ReadUARTData ()

Reads the config string from the device.

Returns

The config string.

Send command to the EPC10

epc10bus	The EPC10 bus
Command	The command

Sets the ADC offset of the DACQ for a single channel

Parameters

AdcChannel	The ADC channel
Offset	The offset for the given channel number

Writes the ADC offset of the DACQ for a single channel to permanent EEProm memory

Parameters

AdcChannel The ADC channel

Sets the parameter of the audio DAC output.

Parameters

Frequency	Frequency(1 - 25000 Hz)
Amplification	Amplification(0 - 0xFFFF)

Sets the DAC Idle value

DacChannel	The DAC channel
IdleValue	The idle value

Sets the offset of a DAC channel.

Parameters

DacChannel	The DAC channel	
Offset	The offset for the given channel number	

Writes the DAC offset of the STG for a single channel to permanent EEProm memory

Parameters

DacChannel	The DAC channel
------------	-----------------

Sets if the DAC Idle value is used after stimulation

Parameters

DacChannel	The DAC channel
Useldle	Use idle value

```
11.43.3.25 SetDigOutState() void SetDigOutState ( uint16_t DigOutState )
```

Writes to the LIH30 digital output

Parameters

DigOutState	The bit mask defining the digital output state

11.43.3.26 SetEEpromPage() void SetEEpromPage (

```
uint32_t EEpromStartAddress,
array< int8_t >^ EEpromData,
LIH30_EPC10_Bus_EnumNet epc10bus )
```

Writes into EEprom on the EPC10 EEPROM

Parameters

EEpromStartAddress	start address of memory area to write to
EEpromData	pointer to internal memory for the supported amount of data
epc10bus	The EPC10 bus

```
11.43.3.27 SetSampleInterval() void SetSampleInterval ( uint32_t SampleInterval )
```

Sets the Sample Interval for the DACQ and Stimulation

Parameters

SampleInterval	between the samples, Sample interval is available from 1 to 4194303
----------------	---

```
11.43.3.28 WriteUARTData() void WriteUARTData (
String^ commandString)
```

Write the command string to the device.

Parameters


commandString	The config string.
---------------	--------------------

11.43.4 Property Documentation

11.43.4.1 StimulusFunction CStimulusFunctionNet^ StimulusFunction [get]

11.44 CMcsBus_AxisParametersNet Class Reference

Inheritance diagram for CMcsBus_AxisParametersNet:

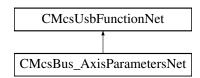

Public Member Functions

- CMcsBus_AxisParametersNet (CMcsUsbNet[^] device)
- ~CMcsBus_AxisParametersNet (void)
- void SetAxisParametersEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, unsigned int parameter)
- void SetAxisParametersEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index, int parameter)
- unsigned int GetAxisParametersUnsignedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)
- int GetAxisParametersSignedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short index)

Additional Inherited Members

11.44.1 Constructor & Destructor Documentation

11.44.2 Member Function Documentation

```
11.44.2.1 GetAxisParametersSignedEeprom() int GetAxisParametersSignedEeprom (
          unsigned char busnumber,
          unsigned char busaddress,
          unsigned char axis,
```

```
11.44.2.2 GetAxisParametersUnsignedEeprom() unsigned int GetAxisParametersUnsignedEeprom ( unsigned char busnumber,
```

```
unsigned char busaddress,
unsigned char axis,
unsigned short index)
```

unsigned short index)

11.44.2.3 SetAxisParametersEeprom() [1/2] void SetAxisParametersEeprom (


```
unsigned char busnumber,
unsigned char busaddress,
unsigned char axis,
unsigned short index,
int parameter)
```

11.44.2.4 SetAxisParametersEeprom() [2/2] void SetAxisParametersEeprom (

```
unsigned char busnumber,
unsigned char busaddress,
unsigned char axis,
unsigned short index,
unsigned int parameter)
```

11.45 CMcsBus_ExtensionNet Class Reference

Inheritance diagram for CMcsBus ExtensionNet:

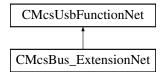

Public Member Functions

- CMcsBus_ExtensionNet (CMcsUsbNet[^] device)
- ~CMcsBus_ExtensionNet (void)
- void SetLEDSwitch (unsigned char busnumber, unsigned char busaddress, unsigned short LEDSwitch)
- unsigned short GetLEDSwitch (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.45.1 Constructor & Destructor Documentation


```
11.45.1.2 \sim CMcsBus_ExtensionNet() \sim CMcsBus_ExtensionNet ( void )
```

11.45.2 Member Function Documentation

11.46 CMcsBus_FYIExtensionNet Class Reference

unsigned char busaddress,
unsigned short LEDSwitch)

Inheritance diagram for CMcsBus_FYIExtensionNet:

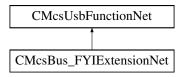

Public Member Functions

- CMcsBus_FYIExtensionNet (CMcsUsbNet[^] device)
- ~CMcsBus_FYIExtensionNet (void)
- void SetValves (unsigned char busnumber, unsigned char busaddress, unsigned int states)
- unsigned int GetValves (unsigned char busnumber, unsigned char busaddress)
- void SetDIO (unsigned char busnumber, unsigned char busaddress, unsigned short io)
- unsigned short GetDIO (unsigned char busnumber, unsigned char busaddress)
- void SetSingleHeater (unsigned char busnumber, unsigned char busaddress, short index, unsigned short power)
- unsigned short GetSingleHeater (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members


11.46.1 Constructor & Destructor Documentation

```
11.46.1.1 CMcsBus_FYIExtensionNet() CMcsBus_FYIExtensionNet ( CMcsUsbNet^ device )
```

```
11.46.1.2 ~CMcsBus_FYIExtensionNet() ~CMcsBus_FYIExtensionNet (
             void )
11.46.2 Member Function Documentation
11.46.2.1 GetDIO() unsigned short GetDIO (
             unsigned char busnumber,
             unsigned char busaddress )
11.46.2.2 GetSingleHeater() unsigned short GetSingleHeater (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.46.2.3 GetValves() unsigned int GetValves (
             unsigned char busnumber,
             unsigned char busaddress )
11.46.2.4 SetDIO() void SetDIO (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short io )
11.46.2.5 SetSingleHeater() void SetSingleHeater (
             unsigned char busnumber,
             unsigned char busaddress,
             short index,
             unsigned short power )
11.46.2.6 SetValves() void SetValves (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned int states )
```

11.47 CMcsBus MotorControlNet Class Reference

Inheritance diagram for CMcsBus_MotorControlNet:

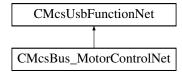

Public Member Functions

- CMcsBus MotorControlNet (CMcsUsbNet[^] device)
- ~CMcsBus MotorControlNet (void)
- void SetMCScalingFactorEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int GetMCScalingFactorEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCScalingFactor (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int factor)
- int GetMCScalingFactor (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short GetMCMaxSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short GetMCMaxSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravelEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int GetMCMaxTravelEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravel (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int GetMCMaxTravel (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- · void SetMCCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCMaxCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- · short GetMCMaxCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRegulatorGainEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short GetMCRegulatorGainEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRegulatorGain (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short gain)
- short GetMCRegulatorGain (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCMaxAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

- void SetMCMaxAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCMaxAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short GetMCStandbyCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short percent)
- short GetMCStandbyCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyTimeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short GetMCStandbyTimeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCStandbyTime (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short t)
- short GetMCStandbyTime (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCBreakCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCBreakCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCBreakCurrent (unsigned char busnumber, unsigned char busnumber, under bus
- short GetMCBreakCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCConfigEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short GetMCConfigEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCConfig (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short config)
- unsigned short GetMCConfig (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short speed)
- unsigned short GetMCSpeedEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short GetMCSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAccelerationEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAcceleration (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReferenceCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCReferenceCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReferenceCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCReferenceCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentModeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet GetMCCurrentModeEeprom (unsigned char busnumber, unsigned char axis)

- void SetMCCurrentMode (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet GetMCCurrentMode (unsigned char busnumber, unsigned char bus
- void SetMCAxisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short revision)
- unsigned short GetMCAxisRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedUnitEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int32_t speedunit)
- int32 t GetMCSpeedUnitEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char axis, bool OnOff_status)
- bool GetMCOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCSpeedShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short speed)
- short GetMCSpeedShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short acceleration)
- unsigned short GetMCAccelerationShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, short current)
- short GetMCCurrentShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCMaxTravelShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int travel)
- int GetMCMaxTravelShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int GetMCCurrentPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCNewPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis, int position)
- int GetMCNewPosition (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- short GetMCCurrentSpeed (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void StartMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCRotation (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char onoff)
- unsigned short GetMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCReference (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned char switch_enable, unsigned char switch_polarity)
- unsigned char GetMCReference (unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned char% switch_port)
- void StopMCMovement (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetMCCurrentModeShortCommand (unsigned char busnumber, unsigned char busaddress, unsigned char axis, RoboCurrentModeEnumNet mode)
- RoboCurrentModeEnumNet GetMCCurrentModeShortCommand (unsigned char busnumber, unsigned char busnumber, unsigned char axis)
- unsigned short GetMCPhase (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- unsigned short GetMCPhaseOffset (unsigned char busnumber, unsigned char busaddress, unsigned char axis)
- void SetSubChannel (unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short subchannel)
- unsigned short GetSubChannel (unsigned char busnumber, unsigned char busaddress, unsigned char axis)

Additional Inherited Members

11.47.1 Constructor & Destructor Documentation

```
11.47.1.1 CMcsBus_MotorControlNet() CMcsBus_MotorControlNet (
             CMcsUsbNet^ device )
11.47.1.2 ~CMcsBus_MotorControlNet() ~CMcsBus_MotorControlNet (
             void )
11.47.2 Member Function Documentation
11.47.2.1 GetMCAcceleration() unsigned short GetMCAcceleration (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.2 GetMCAccelerationEeprom() unsigned short GetMCAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.3 GetMCAccelerationShortCommand() unsigned short GetMCAccelerationShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.4 GetMCAxisRevisionEeprom() unsigned short GetMCAxisRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.5 GetMCBreakCurrent() short GetMCBreakCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.6 GetMCBreakCurrentEeprom() short GetMCBreakCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.7 GetMCConfig() unsigned short GetMCConfig (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
\textbf{11.47.2.8} \quad \textbf{GetMCConfigEeprom()} \quad \texttt{unsigned short GetMCConfigEeprom ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.9 GetMCCurrent() short GetMCCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.10 GetMCCurrentEeprom() short GetMCCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.11 GetMCCurrentMode() RoboCurrentModeEnumNet GetMCCurrentMode (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.12 GetMCCurrentModeEeprom() RoboCurrentModeEnumNet GetMCCurrentModeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.13 GetMCCurrentModeShortCommand() RoboCurrentModeEnumNet GetMCCurrentModeShort←
Command (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.14 GetMCCurrentPosition() int GetMCCurrentPosition (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.15 GetMCCurrentShortCommand() short GetMCCurrentShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.16 GetMCCurrentSpeed() short GetMCCurrentSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
\textbf{11.47.2.17} \quad \textbf{GetMCMaxAcceleration()} \quad \textbf{unsigned short GetMCMaxAcceleration ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.18 GetMCMaxAccelerationEeprom() unsigned short GetMCMaxAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.19 GetMCMaxCurrent() short GetMCMaxCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.20 GetMCMaxCurrentEeprom() short GetMCMaxCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
{\bf 11.47.2.21} \quad {\bf GetMCMaxSpeed()} \quad {\tt unsigned \ short \ GetMCMaxSpeed} \quad {\tt (}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.22 GetMCMaxSpeedEeprom() unsigned short GetMCMaxSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.23 GetMCMaxTravel() int GetMCMaxTravel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.24 GetMCMaxTravelEeprom() int GetMCMaxTravelEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.25 GetMCMaxTravelShortCommand() int GetMCMaxTravelShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.26 GetMCMovement() unsigned short GetMCMovement (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
\textbf{11.47.2.27} \quad \textbf{GetMCNewPosition()} \quad \texttt{int GetMCNewPosition ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.28 GetMCOutputOnOff() bool GetMCOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.29 GetMCPhase() unsigned short GetMCPhase (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.30 GetMCPhaseOffset() unsigned short GetMCPhaseOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.31 GetMCReference() unsigned char GetMCReference (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned char% switch_port)
11.47.2.32 GetMCReferenceCurrent() short GetMCReferenceCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.33 GetMCReferenceCurrentEeprom() short GetMCReferenceCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.34 GetMCRegulatorGain() short GetMCRegulatorGain (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.35 GetMCRegulatorGainEeprom() short GetMCRegulatorGainEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.36 GetMCScalingFactor() int GetMCScalingFactor (
             unsigned char busnumber,
             unsigned char busaddress,
            unsigned char axis )
11.47.2.37 GetMCScalingFactorEeprom() int GetMCScalingFactorEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.38 GetMCSpeed() short GetMCSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.39 GetMCSpeedEeprom() unsigned short GetMCSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.40 GetMCSpeedShortCommand() short GetMCSpeedShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.41 GetMCSpeedUnitEeprom() int32_t GetMCSpeedUnitEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.42 GetMCStandbyCurrent() short GetMCStandbyCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.43 GetMCStandbyCurrentEeprom() short GetMCStandbyCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.44 GetMCStandbyTime() short GetMCStandbyTime (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.45 GetMCStandbyTimeEeprom() short GetMCStandbyTimeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
11.47.2.46 GetSubChannel() unsigned short GetSubChannel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis )
```

```
11.47.2.47 SetMCAcceleration() void SetMCAcceleration (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.47.2.48 SetMCAccelerationEeprom() void SetMCAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.47.2.49 SetMCAccelerationShortCommand() void SetMCAccelerationShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.47.2.50 SetMCAxisRevisionEeprom() void SetMCAxisRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short revision )
11.47.2.51 SetMCBreakCurrent() void SetMCBreakCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.52 SetMCBreakCurrentEeprom() void SetMCBreakCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
```

```
11.47.2.53 SetMCConfig() void SetMCConfig (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short config )
11.47.2.54 SetMCConfigEeprom() void SetMCConfigEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short config )
11.47.2.55 SetMCCurrent() void SetMCCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.56 SetMCCurrentEeprom() void SetMCCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.57 SetMCCurrentMode() void SetMCCurrentMode (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             RoboCurrentModeEnumNet mode )
11.47.2.58 SetMCCurrentModeEeprom() void SetMCCurrentModeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             RoboCurrentModeEnumNet mode )
```

```
11.47.2.59 SetMCCurrentModeShortCommand() void SetMCCurrentModeShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             RoboCurrentModeEnumNet mode )
11.47.2.60 SetMCCurrentPosition() void SetMCCurrentPosition (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int position )
11.47.2.61 SetMCCurrentShortCommand() void SetMCCurrentShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
\textbf{11.47.2.62} \quad \textbf{SetMCMaxAcceleration()} \quad \texttt{void SetMCMaxAcceleration ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.47.2.63 SetMCMaxAccelerationEeprom() void SetMCMaxAccelerationEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short acceleration )
11.47.2.64 SetMCMaxCurrent() void SetMCMaxCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
```

```
11.47.2.65 SetMCMaxCurrentEeprom() void SetMCMaxCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.66 SetMCMaxSpeed() void SetMCMaxSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
11.47.2.67 SetMCMaxSpeedEeprom() void SetMCMaxSpeedEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
11.47.2.68 SetMCMaxTravel() void SetMCMaxTravel (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
11.47.2.69 SetMCMaxTravelEeprom() void SetMCMaxTravelEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
11.47.2.70 SetMCMaxTravelShortCommand() void SetMCMaxTravelShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int travel )
```


```
11.47.2.71 SetMCNewPosition() void SetMCNewPosition (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int position )
11.47.2.72 SetMCOutputOnOff() void SetMCOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             bool OnOff_status )
11.47.2.73 SetMCReference() void SetMCReference (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned char switch_enable,
             unsigned char switch_polarity)
11.47.2.74 SetMCReferenceCurrent() void SetMCReferenceCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.75 SetMCReferenceCurrentEeprom() void SetMCReferenceCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short current )
11.47.2.76 SetMCRegulatorGain() void SetMCRegulatorGain (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short gain )
```

```
11.47.2.77 SetMCRegulatorGainEeprom() void SetMCRegulatorGainEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short gain )
11.47.2.78 SetMCRotation() void SetMCRotation (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned char onoff )
11.47.2.79 SetMCScalingFactor() void SetMCScalingFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int factor )
11.47.2.80 SetMCScalingFactorEeprom() void SetMCScalingFactorEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int factor )
11.47.2.81 SetMCSpeed() void SetMCSpeed (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short speed )
\textbf{11.47.2.82} \quad \textbf{SetMCSpeedEeprom()} \quad \texttt{void SetMCSpeedEeprom ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short speed )
```

```
11.47.2.83 SetMCSpeedShortCommand() void SetMCSpeedShortCommand (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short speed )
11.47.2.84 SetMCSpeedUnitEeprom() void SetMCSpeedUnitEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int32_t speedunit )
11.47.2.85 SetMCStandbyCurrent() void SetMCStandbyCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short percent )
11.47.2.86 SetMCStandbyCurrentEeprom() void SetMCStandbyCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short percent )
11.47.2.87 SetMCStandbyTime() void SetMCStandbyTime (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short t )
11.47.2.88 SetMCStandbyTimeEeprom() void SetMCStandbyTimeEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short t)
```

11.48 CMcsBus SensorNet Class Reference

Inheritance diagram for CMcsBus_SensorNet:

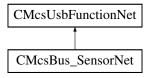

Public Member Functions

- CMcsBus SensorNet (CMcsUsbNet[^] device)
- ~CMcsBus_SensorNet (void)
- void SetMinimalThreshold (unsigned char busnumber, unsigned char busaddress, unsigned short threshold)
- unsigned short GetMinimalThreshold (unsigned char busnumber, unsigned char busaddress)
- void SetDetectionThreshold (unsigned char busnumber, unsigned char busnumber, unsigned short threshold)
- unsigned short GetDetectionThreshold (unsigned char busnumber, unsigned char busaddress)
- void SetLatency (unsigned char busnumber, unsigned char busaddress, unsigned short latency)
- unsigned short GetLatency (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetBubbleStatus (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetLatencyCounter (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetDetectorValue (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ GetPressure (unsigned char busnumber, unsigned char busaddress, int n)
- int GetPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRegulatorOnOff (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned char onoff)
- unsigned char GetRegulatorOnOff (unsigned char busnumber, unsigned char busaddress, unsigned short index)

- void SetSollPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index, int pressure)
- int GetSollPressure (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRegulatorFactor (unsigned char busnumber, unsigned char busaddress, unsigned short index, int factor)
- int GetRegulatorFactor (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetPressureOffset (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- array< unsigned short > ^ GetPressureOffset (unsigned char busnumber, unsigned char busaddress)
- int GetPressureOffset (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- unsigned int GetRegulatorStatus (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetRotatePump (unsigned char busnumber, unsigned char busaddress, unsigned short index, short speed)
- short GetRotatePump (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetMovePump (unsigned char busnumber, unsigned char busaddress, unsigned short index, unsigned short speed, int position)
- void SetRegulationTimeouts (unsigned char busnumber, unsigned char busaddress, unsigned short Max
 SpeedWait, unsigned short MaxSignChange)
- array< int > ^ Get4ADC (unsigned char busnumber, unsigned char busaddress)
- array< int > ^ Get4ADCAverage (unsigned char busnumber, unsigned char busaddress)
- void Set4DAC (unsigned char busnumber, unsigned char busaddress, array< unsigned short >^ dac)
- array< unsigned short > ^ Get4DAC (unsigned char busnumber, unsigned char busaddress)
- void Set4ADCMode (unsigned char busnumber, unsigned char busaddress, PatchServAdcModeEnumNet mode)
- PatchServAdcModeEnumNet Get4ADCMode (unsigned char busnumber, unsigned char busaddress)
- void Set4ADCCatchampAverageShift (unsigned char busnumber, unsigned char busaddress, unsigned int shift)
- unsigned int Get4ADCCatchampAverageShift (unsigned char busnumber, unsigned char busaddress)
- array < unsigned short > ^ Get2AnalogInput (unsigned char busnumber, unsigned char busaddress)
- unsigned short Get2DigitalInput (unsigned char busnumber, unsigned char busaddress)
- array< unsigned short > ^ GetADCs (unsigned char busnumber, unsigned char busaddress, int n)
- array< unsigned short > $^{\wedge}$ GetADCsLoop (unsigned char busnumber, unsigned char busaddress, int n)
- void SetPiezoState (unsigned char busnumber, unsigned char busaddress, int state)
- void SetDACs (unsigned char busnumber, unsigned char busaddress, unsigned short index, array< unsigned short >^ dac_times_voltages)
- array< unsigned short > ^ GetDACs (unsigned char busnumber, unsigned char busaddress, unsigned short index)
- void SetSamplePeriode (unsigned char busnumber, unsigned char busaddress, unsigned short periode)
- unsigned short GetSamplePeriode (unsigned char busnumber, unsigned char busaddress)
- void StartSync (unsigned char busnumber, unsigned char busaddress)
- unsigned short GetSyncState (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacAmplitude (unsigned char busnumber, unsigned char busaddress, unsigned short dacAmplitude)
- unsigned short CatchAmpGetDacAmplitude (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacOffset (unsigned char busnumber, unsigned char busaddress, short dacOffset)
- short CatchAmpGetDacOffset (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcMean (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcValue (unsigned char busnumber, unsigned char busaddress)

- int CatchAmpGetAdcValueH (unsigned char busnumber, unsigned char busaddress)
- int CatchAmpGetAdcValueL (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetPwmEnable (unsigned char busnumber, unsigned char busaddress, bool pwmEnable)
- bool CatchAmpGetPwmEnable (unsigned char busnumber, unsigned char busaddress)
- void CatchAmpSetDacEnable (unsigned char busnumber, unsigned char busaddress, bool dacEnable)
- bool CatchAmpGetDacEnable (unsigned char busnumber, unsigned char busaddress)
- int TactSwitchGetState (unsigned char busnumber, unsigned char busaddress)
- void TactSwitchSetDisplay (unsigned char busnumber, unsigned char busaddress, int Melody)

Additional Inherited Members

11.48.1 Constructor & Destructor Documentation

11.48.2 Member Function Documentation

```
11.48.2.1 CatchAmpGetAdcMean() int CatchAmpGetAdcMean (
unsigned char busnumber,
unsigned char busaddress)
```

```
11.48.2.2 CatchAmpGetAdcValue() int CatchAmpGetAdcValue (
unsigned char busnumber,
unsigned char busaddress)
```

```
11.48.2.4 CatchAmpGetAdcValueL() int CatchAmpGetAdcValueL (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.5 CatchAmpGetDacAmplitude() unsigned short CatchAmpGetDacAmplitude (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.6 CatchAmpGetDacEnable() bool CatchAmpGetDacEnable (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.7 CatchAmpGetDacOffset() short CatchAmpGetDacOffset (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.8 CatchAmpGetPwmEnable() bool CatchAmpGetPwmEnable (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.9 CatchAmpSetDacAmplitude() void CatchAmpSetDacAmplitude (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short dacAmplitude )
11.48.2.10 CatchAmpSetDacEnable() void CatchAmpSetDacEnable (
             unsigned char busnumber,
             unsigned char busaddress,
             bool dacEnable )
11.48.2.11 CatchAmpSetDacOffset() void CatchAmpSetDacOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             short dacOffset )
```

```
11.48.2.12 CatchAmpSetPwmEnable() void CatchAmpSetPwmEnable (
             unsigned char busnumber,
             unsigned char busaddress,
             bool pwmEnable )
11.48.2.13 Get2AnalogInput() array<unsigned short> ^ Get2AnalogInput (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.14 Get2DigitalInput() unsigned short Get2DigitalInput (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.15 Get4ADC() array<int> ^ Get4ADC (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.16 Get4ADCAverage() array<int> ^ Get4ADCAverage (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.17 Get4ADCCatchampAverageShift() unsigned int Get4ADCCatchampAverageShift (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.18 Get4ADCMode() PatchServAdcModeEnumNet Get4ADCMode (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.19 Get4DAC() array<unsigned short> ^{\land} Get4DAC (
             unsigned char busnumber,
             unsigned char busaddress )
```

```
11.48.2.20 GetADCs() array<unsigned short> ^ GetADCs (
             unsigned char busnumber,
             unsigned char busaddress,
             int n)
11.48.2.21 GetADCsLoop() array<unsigned short> ^ GetADCsLoop (
            unsigned char busnumber,
             unsigned char busaddress,
             int n)
11.48.2.22 GetBubbleStatus() unsigned short GetBubbleStatus (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.23 GetDACs() array<unsigned short> ^ GetDACs (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
11.48.2.24 GetDetectionThreshold() unsigned short GetDetectionThreshold (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.25 GetDetectorValue() unsigned short GetDetectorValue (
             unsigned char busnumber,
            unsigned char busaddress )
11.48.2.26 GetLatency() unsigned short GetLatency (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.27 GetLatencyCounter() unsigned short GetLatencyCounter (
             unsigned char busnumber,
             unsigned char busaddress )
```

```
11.48.2.28 GetMinimalThreshold() unsigned short GetMinimalThreshold (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.29 GetMovePump() void GetMovePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned short% speed,
             [System::Runtime::InteropServices::Out] int% position )
11.48.2.30 GetPiezoState() void GetPiezoState (
             unsigned char busnumber,
             unsigned char busaddress,
             [System::Runtime::InteropServices::Out] int% state,
             [System::Runtime::InteropServices::Out] int% reason )
11.48.2.31 GetPressure() [1/2] array<int> ^ GetPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             int n)
11.48.2.32 GetPressure() [2/2] int GetPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.48.2.33 GetPressureOffset() [1/2] array<unsigned short> ^ GetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.34 GetPressureOffset() [2/2] int GetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
```

```
11.48.2.35 GetRegulationTimeouts() void GetRegulationTimeouts (
             unsigned char busnumber,
             unsigned char busaddress,
             [System::Runtime::InteropServices::Out] unsigned short% MaxSpeedWait,
             [System::Runtime::InteropServices::Out] unsigned short% MaxSignChange )
11.48.2.36 GetRegulatorFactor() int GetRegulatorFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.48.2.37 GetRegulatorOnOff() unsigned char GetRegulatorOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.48.2.38 GetRegulatorStatus() unsigned int GetRegulatorStatus (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index)
11.48.2.39 GetRotatePump() short GetRotatePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
{\bf 11.48.2.40} \quad {\bf GetSamplePeriode()} \quad {\tt unsigned short GetSamplePeriode ()}
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.41 GetSollPressure() int GetSollPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
```


```
11.48.2.42 GetSyncState() unsigned short GetSyncState (
              unsigned char busnumber,
              unsigned char busaddress )
\textbf{11.48.2.43} \quad \textbf{Set4ADCCatchampAverageShift()} \quad \texttt{void Set4ADCCatchampAverageShift ()}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned int shift )
11.48.2.44 Set4ADCMode() void Set4ADCMode (
              unsigned char busnumber,
              unsigned char busaddress,
              PatchServAdcModeEnumNet mode )
11.48.2.45 Set4DAC() void Set4DAC (
              unsigned char busnumber,
              unsigned char busaddress,
              array < unsigned short >^{\wedge} dac)
11.48.2.46 SetDACs() void SetDACs (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short index,
              array< unsigned short >^{\wedge} dac\_times\_voltages )
\textbf{11.48.2.47} \quad \textbf{SetDetectionThreshold()} \quad \texttt{void SetDetectionThreshold ()}
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short threshold )
11.48.2.48 SetLatency() void SetLatency (
              unsigned char busnumber,
              unsigned char busaddress,
              unsigned short latency )
```

```
11.48.2.49 SetMinimalThreshold() void SetMinimalThreshold (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short threshold )
11.48.2.50 SetMovePump() void SetMovePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             unsigned short speed,
             int position )
11.48.2.51 SetPiezoState() void SetPiezoState (
             unsigned char busnumber,
             unsigned char busaddress,
             int state )
11.48.2.52 SetPressureOffset() void SetPressureOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index )
11.48.2.53 SetRegulationTimeouts() void SetRegulationTimeouts (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short MaxSpeedWait,
             unsigned short MaxSignChange )
11.48.2.54 SetRegulatorFactor() void SetRegulatorFactor (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             int factor )
11.48.2.55 SetRegulatorOnOff() void SetRegulatorOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             unsigned char onoff )
```

```
11.48.2.56 SetRotatePump() void SetRotatePump (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             short speed )
11.48.2.57 SetSamplePeriode() void SetSamplePeriode (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short periode )
11.48.2.58 SetSollPressure() void SetSollPressure (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short index,
             int pressure )
11.48.2.59 StartSync() void StartSync (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.60 TactSwitchGetState() int TactSwitchGetState (
             unsigned char busnumber,
             unsigned char busaddress )
11.48.2.61 TactSwitchSetDisplay() void TactSwitchSetDisplay (
             unsigned char busnumber,
             unsigned char busaddress,
             int Melody )
```

11.49 CMcsBus_TempSensorNet Class Reference

Inheritance diagram for CMcsBus_TempSensorNet:

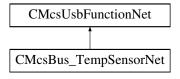

Public Member Functions

- CMcsBus_TempSensorNet (CMcsUsbNet[^] device)
- ~CMcsBus_TempSensorNet (void)
- short GetTemperatur (unsigned char busnumber, unsigned char busaddress)
- short GetTemperatur (unsigned char busnumber, unsigned char busaddress, short index)
- void SetNanoVoltsPerKelvin (unsigned char busnumber, unsigned char busaddress, int nanovoltsperkelvin)
- int GetNanoVoltsPerKelvin (unsigned char busnumber, unsigned char busaddress)
- short GetThermoVoltage (unsigned char busnumber, unsigned char busaddress, short index)
- short GetThermoTemp (unsigned char busnumber, unsigned char busaddress, short index)
- · void SetThermoOffset (unsigned char busnumber, unsigned char busaddress, short index, short offset)
- short GetThermoOffset (unsigned char busnumber, unsigned char busaddress, short index)

Additional Inherited Members


11.49.1 Constructor & Destructor Documentation

short index)

```
11.49.2.4 GetThermoOffset() short GetThermoOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.49.2.5 GetThermoTemp() short GetThermoTemp (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.49.2.6 GetThermoVoltage() short GetThermoVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             short index )
11.49.2.7 SetNanoVoltsPerKelvin() void SetNanoVoltsPerKelvin (
             unsigned char busnumber,
             unsigned char busaddress,
             int nanovoltsperkelvin )
11.49.2.8 SetThermoOffset() void SetThermoOffset (
             unsigned char busnumber,
             unsigned char busaddress,
             short index,
             short offset )
```

11.50 CMcsBus_VoltageModeNet Class Reference

Inheritance diagram for CMcsBus_VoltageModeNet:

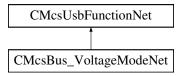

Public Member Functions

- CMcsBus VoltageModeNet (CMcsUsbNet[^] device)
- ~CMcsBus_VoltageModeNet (void)
- void SetVMMaxPositiveCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxPositiveCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxPositiveCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxNegativeCurrentEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short current)
- short GetVMMaxNegativeCurrent (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxPositiveVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxPositiveVoltage (unsigned char busnumber, unsigned char
- short GetVMMaxPositiveVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxNegativeVoltageEeprom (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMMaxNegativeVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMMaxNegativeVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char channel, unsigned short status)
- unsigned short GetVMOutputOnOff (unsigned char busnumber, unsigned char busaddress, unsigned char channel)
- void SetVMVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel, short voltage)
- short GetVMVoltage (unsigned char busnumber, unsigned char busaddress, unsigned char channel)

Additional Inherited Members

11.50.1 Constructor & Destructor Documentation



```
11.50.1.2 ~ CMcsBus_VoltageModeNet() ~ CMcsBus_VoltageModeNet (
             void )
11.50.2 Member Function Documentation
11.50.2.1 GetVMMaxNegativeCurrent() short GetVMMaxNegativeCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.2 GetVMMaxNegativeCurrentEeprom() short GetVMMaxNegativeCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.3 GetVMMaxNegativeVoltage() short GetVMMaxNegativeVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char {\it channel} )
11.50.2.4 GetVMMaxNegativeVoltageEeprom() short GetVMMaxNegativeVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
\textbf{11.50.2.5} \quad \textbf{GetVMMaxPositiveCurrent()} \quad \texttt{short GetVMMaxPositiveCurrent ()}
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.6 GetVMMaxPositiveCurrentEeprom() short GetVMMaxPositiveCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
```

```
11.50.2.7 GetVMMaxPositiveVoltage() short GetVMMaxPositiveVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.8 GetVMMaxPositiveVoltageEeprom() short GetVMMaxPositiveVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.9 GetVMOutputOnOff() unsigned short GetVMOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.10 GetVMVoltage() short GetVMVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel )
11.50.2.11 SetVMMaxNegativeCurrent() void SetVMMaxNegativeCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.50.2.12 SetVMMaxNegativeCurrentEeprom() void SetVMMaxNegativeCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.50.2.13 SetVMMaxNegativeVoltage() void SetVMMaxNegativeVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
```

```
11.50.2.14 SetVMMaxNegativeVoltageEeprom() void SetVMMaxNegativeVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
11.50.2.15 SetVMMaxPositiveCurrent() void SetVMMaxPositiveCurrent (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.50.2.16 SetVMMaxPositiveCurrentEeprom() void SetVMMaxPositiveCurrentEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short current )
11.50.2.17 SetVMMaxPositiveVoltage() void SetVMMaxPositiveVoltage (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
11.50.2.18 SetVMMaxPositiveVoltageEeprom() void SetVMMaxPositiveVoltageEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
11.50.2.19 SetVMOutputOnOff() void SetVMOutputOnOff (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             unsigned short status )
 \textbf{11.50.2.20} \quad \textbf{SetVMVoltage()} \quad \texttt{void SetVMVoltage ()} 
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char channel,
             short voltage )
```

11.51 CMcsBusNet Class Reference

Inheritance diagram for CMcsBusNet:

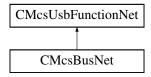

Public Member Functions

- CMcsBusNet (CMcsUsbNet[^] device)
- virtual ∼CMcsBusNet (void)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned short value)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, short value)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, unsigned int value)
- void SetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, int value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned short% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]short% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]unsigned int% value)
- void GetCommand (unsigned char command, unsigned char busnumber, unsigned char busaddress, unsigned char axis, [System::Runtime::InteropServices::Out]int% value)
- void SetBusAddressEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short GetBusAddressEeprom (unsigned char busnumber, unsigned char busaddress)
- void SetBusAddress (unsigned char busnumber, unsigned char busaddress, unsigned short newaddress)
- unsigned short GetBusAddress (unsigned char busnumber, unsigned char busaddress)
- void CMcsBusNet::SetModeEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short CMcsBusNet::GetModeEeprom (unsigned char busnumber, unsigned char busaddress)
- void CMcsBusNet::SetMode (unsigned char busnumber, unsigned char busaddress, unsigned short mode)
- unsigned short CMcsBusNet::GetMode (unsigned char busnumber, unsigned char busaddress)
- void SetHWRevisionEeprom (unsigned char busnumber, unsigned char busaddress, unsigned short revision)
- unsigned short GetHWRevisionEeprom (unsigned char busnumber, unsigned char busaddress)

Additional Inherited Members

11.51.1 Constructor & Destructor Documentation

```
11.51.1.1 CMcsBusNet() CMcsBusNet (
CMcsUsbNet^ device)
```

```
11.51.1.2 ~CMcsBusNet() virtual ~CMcsBusNet (
            void ) [virtual]
11.51.2 Member Function Documentation
11.51.2.1 CMcsBusNet::GetMode() unsigned short CMcsBusNet::GetMode (
            unsigned char busnumber,
            unsigned char busaddress )
11.51.2.2 CMcsBusNet::GetModeEeprom() unsigned short CMcsBusNet::GetModeEeprom (
            unsigned char busnumber,
            unsigned char busaddress )
11.51.2.3 CMcsBusNet::SetMode() void CMcsBusNet::SetMode (
            unsigned char busnumber,
            unsigned char busaddress,
            unsigned short mode )
11.51.2.4 CMcsBusNet::SetModeEeprom() void CMcsBusNet::SetModeEeprom (
            unsigned char busnumber,
            unsigned char busaddress,
            unsigned short mode )
11.51.2.5 GetBusAddress() unsigned short GetBusAddress (
            unsigned char busnumber,
            unsigned char busaddress )
11.51.2.6 GetBusAddressEeprom() unsigned short GetBusAddressEeprom (
            unsigned char busnumber,
            unsigned char busaddress )
```


```
11.51.2.7 GetCommand() [1/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] int% value )
11.51.2.8 GetCommand() [2/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] short% value )
11.51.2.9 GetCommand() [3/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned int% value )
11.51.2.10 GetCommand() [4/4] void GetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             [System::Runtime::InteropServices::Out] unsigned short% value )
11.51.2.11 GetHWRevisionEeprom() unsigned short GetHWRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress )
11.51.2.12 SetBusAddress() void SetBusAddress (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short newaddress )
```

```
11.51.2.13 SetBusAddressEeprom() void SetBusAddressEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short newaddress )
11.51.2.14 SetCommand() [1/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             int value )
11.51.2.15 SetCommand() [2/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             short value )
11.51.2.16 SetCommand() [3/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned int value )
11.51.2.17 SetCommand() [4/4] void SetCommand (
             unsigned char command,
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned char axis,
             unsigned short value )
11.51.2.18 SetHWRevisionEeprom() void SetHWRevisionEeprom (
             unsigned char busnumber,
             unsigned char busaddress,
             unsigned short revision )
```

11.52 CMcsUsbDacqNet Class Reference

Base class for data acquisition devices.

Inheritance diagram for CMcsUsbDacqNet:

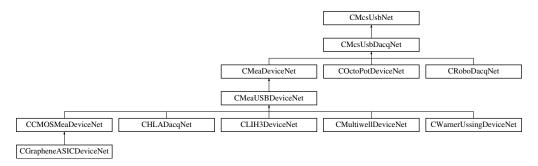

Classes

· class CHWInfo

Class to provide hardware information about the device.

Public Member Functions

- CMcsUsbDacqNet ()
- ∼CMcsUsbDacqNet ()
- uint32_t GetErrorMessage ([System::Runtime::InteropServices::Out]String[^]% errorString, [System::←
 Runtime::InteropServices::Out]int% info)
- virtual uint32_t GetVoltageRangeIndex (unsigned int virtualDevice)
- virtual void SetVoltageRangeByIndex (int32 t voltageRangeIndex, unsigned int virtualDevice)

Sets the voltage range on devices which support multiple voltage ranges.

virtual void SetVoltageRangeInMicroVolt (int32_t voltageRange, unsigned int virtualDevice)

Sets the voltage range on devices which support multiple voltage ranges.

virtual int32_t GetVoltageRangeInMicroVolt (unsigned int virtualDevice)

Gets the currently selected voltage range on devices which support multiple voltage ranges.

virtual int32_t GetVoltageRangeInMilliVolt ()

Gets the currently selected voltage range on devices which support multiple voltage ranges.

virtual void SetDataMode (DataModeEnumNet dataMode, unsigned int virtualDevice)

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

virtual DataModeEnumNet GetDataMode (unsigned int virtualDevice)

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MEA2100DigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

• void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MultiwellDigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

• void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, CMOSMEA5000DigitalSourceEnumNet source, int bitnumber offset)

Sets the function/source of an digital output bit.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, W2100DigitalSourceEnumNet source, int bitnumber_offset) Sets the function/source of an digital output bit.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, SCUDigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, MEA2100_256DigitalSourceEnumNet source, int bitnumber_offset)

Sets the function/source of an digital output bit.

template<typename digitalsourceenum >
 void SetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, DigitalSource< digitalsourceenum >^ source, int bitnumber offset)

Sets the function/source of an digital output bit.

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::
 — InteropServices::Out]MultiwellDigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::
 — InteropServices::Out]W2100DigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber offset)

Gets the function/source of an digital output bit.

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::
 InteropServices::Out]SCUDigitalSourceEnumNet% source, [System::Runtime::InteropServices::Out]int% bitnumber_offset)

Gets the function/source of an digital output bit.

Gets the function/source of an digital output bit.

template<typename digitalsourceenum >

void GetDigitalSource (DigitalTargetEnumNet digitaltarget, int32_t NrChannel, [System::Runtime::Interop⇔ Services::Out]DigitalSource< digitalsourceenum >^% source, [System::Runtime::InteropServices::Out]int% bitnumber offset)

Gets the function/source of an digital output bit.

virtual AdapterTypeEnumNet GetAdapterType ()

Gets the adapter which is connected to the MEA2100 device.

virtual MeaLayoutEnumNet GetMeaLayout ()

Gets the MEA layout which is connected to the MEA2100 device.

virtual uint32_t GetAdcDataFormat (uint32_t virtualDevice)

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

- virtual uint32_t GetResolutionPerDigit (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% res, [System::Runtime::InteropServices::Out] int% resUnit)
- virtual uint32_t GetHardwareMinRange (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)

- virtual uint32_t GetHardwareMaxRange (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% r, [System::Runtime::InteropServices::Out] int% rUnit)
- virtual uint32_t GetNumberOfDataBits (uint32_t virtualDevice, DacqGroupChannelEnumNet group, [System::Runtime::InteropServices::Out] int% numberOfBits)

Get the real number of data bits.

virtual void SetSamplerate (int32_t rate, unsigned int oversample, unsigned int virtualDevice)

Sets the sampling frequency of the device.

virtual int32_t GetSamplerate (unsigned int virtualDevice)

Gets the sampling frequency of the device.

virtual uint32_t GetMaxSamplingFrequency (int virtualDevice)

Gets the maximal sampling frequency of the device.

virtual uint32_t GetMinSamplingFrequencyStepsize ()

Gets the minimal sampling frequency step size increment value of the device.

• virtual int32_t GetChannelsInBlock (unsigned int virtualDevice)

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

- virtual void GetChannelLayout ([System::Runtime::InteropServices::Out]int% AnalogChannels, [System::Runtime::InteropServices::Out]int% DigitalChannels, [System::Runtime::InteropServices::Out]int% ChecksumChannels, [System::Runtime::InteropServices::Out]int% TimestampChannels, [System::← Runtime::InteropServices::Out]int% ChannelsInBlock, unsigned int virtualDevice)
- virtual void SendStartDacq ()

Start sampling.

virtual void SendStartDacq (int VirtualDacqMap)

Start sampling

virtual void SendStartStgAndDacq (uint32_t trigger_map, int VirtualDacqMap)

Start sampling together with the STG.

virtual void SendStopDacq ()

Stop sampling.

virtual void SendStopDacq (int VirtualDacqMap)

Stop sampling.

Parameters

VirtualDacqMap

virtual void SendStopStgAndDacq (uint32_t trigger_map, int VirtualDacqMap)

Stop sampling together with the STG.

virtual void SendStopStgAndDacqWithOptions (uint32_t trigger_map, int VirtualDacqMap, int options)

Stop sampling together with the STG and options.

virtual void StartLoop ()

Start the data acquisition thread.

virtual void StartLoop (int32_t timeout)

Start the data acquisition thread.

virtual void StartLoop (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_← t packetsInUrb)

Start the data acquisition thread.

 virtual void StartLoop (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_← t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread.

virtual void StopLoop ()

- · virtual void ClearBuffers ()
- virtual void StartDacq ()

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout)

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_← t packetsInUrb)

Start the data acquisition thread and sampling.

virtual void StartDacq (int32_t timeout, int32_t numSubmittedUsbBuffers, int32_t numUsbBuffers, int32_←
t packetsInUrb, uint32_t virtualDevice)

Start the data acquisition thread and sampling.

virtual void StopDacq ()

Stop the data acquisition thread and sampling.

virtual void StopDacq (uint32_t virtualDevice)

Stop the data acquisition thread and sampling.

- virtual uint32_t SetPoti (uint32_t channel, uint32_t value, bool write_nvram)
- virtual uint32 t GetPoti (uint32 t channel, [System::Runtime::InteropServices::Out]uint32 t% value)
- virtual CFilterPropertyNet ^ GetFilterProperty (DacqGroupChannelEnumNet GroupID, unsigned int index)
- virtual array
 CFilterPropertyNet[^]> [^] CMcsUsbDacqNet::GetFilterProperties (DacqGroupChannelEnumNet GroupID)
- int GetChannelDataFillSize ()
- virtual void SetSelectedChannels (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void SetSelectedChannels (int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedChannels (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

- virtual void SetSelectedChannels (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedData (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void SetSelectedData (int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedData (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

- virtual void SetSelectedData (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampl
- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read← FramesDict... with handle = 0 to read the data.

- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, int nChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSizeNet sampleDstSize)
- virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, array< bool >[^] selected←
 Channels, int queuesize, int threshold, SampleSizeNet samplesize)

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read←FramesDict... with handle = 0 to read the data.

virtual int AddSelectedChannelsQueue (int nByteOffset, int nChannelOffset, array< bool >[^] selected←
 Channels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleDstSizeNet

 virtual void SetSelectedChannelsQueue (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

- virtual void SetSelectedChannelsQueue (int nChannels, int queuesize, int threshold, SampleSizeNet samplesize, SampleDstSizeNet sampleDstSize, int ChannelsInBlock)
- virtual void SetSelectedChannelsQueue (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet samplesize, int ChannelsInBlock)

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

- virtual void SetSelectedChannelsQueue (array< bool >^ selectedChannels, int queuesize, int threshold, SampleSizeNet sampleDstSizeNet sampleD
- virtual uint32 t ChannelBlock AvailFrames (int handle)

Get the number of sample frames already available in the FIFO.

- virtual uint32_t ChannelBlock_AvailFrames (int handle, int queue)

Read data from a FIFO queue in uint16_t data format

• virtual void ChannelBlock_ReadFramesUI16 (int handle, array< uint16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format

virtual array< int16_t > ^ ChannelBlock_ReadFramesI16 (int handle, int frames, [System::Runtime::
 —
 InteropServices::Out]int % frames ret)

Read data from a FIFO queue in int16_t data format

virtual void ChannelBlock_ReadFramesI16 (int handle, array< int16_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format

Read data from a FIFO queue in uint32_t data format

• virtual void ChannelBlock_ReadFramesUI32 (int handle, array< uint32_t >^ buffer, int frames_pos, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

• virtual array< int32_t > ^ ChannelBlock_ReadFramesI32 (int handle, int frames, [System::Runtime:: InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint32_t data format

virtual void ChannelBlock_ReadFramesl32 (int handle, array< int32_t >^ buffer, int frames_pos, int frames,
 [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format

• virtual array< array< uint16_t $>^{^{\wedge}}>$ ^ ChannelBlock_ReadAsFrameArrayUI16 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< uint16_t >^> ^ ChannelBlock_ReadAsFrameArrayUI16 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

virtual array< array< int16_t >^> ^ ChannelBlock_ReadAsFrameArrayI16 (int handle, int frames, [System← ::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

virtual array< array< int16_t >^> ^ ChannelBlock_ReadAsFrameArrayI16 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

virtual array< array< uint32_t >^> ^ ChannelBlock_ReadAsFrameArrayUl32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< uint32_t >^> ^ ChannelBlock_ReadAsFrameArrayUI32 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

virtual array < array < int32_t >^ > ^ ChannelBlock_ReadAsFrameArrayl32 (int handle, int frames, [System ← ::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

• virtual array< array< int32_t >^> ^ ChannelBlock_ReadAsFrameArrayl32 (int handle, int queue, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue as array of uint16_t data frame arrays

 virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ ChannelBlock_ReadFramesDictUI16 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ ChannelBlock_ReadFramesDict116
 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ ChannelBlock_ReadFramesDictUI32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

virtual System::Collections::Generic::Dictionary< int, array< int32_t >^> ^ ChannelBlock_ReadFramesDictl32 (int handle, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint16_t >^> ^ GetGroupChannelDataUI16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames ret)

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

 virtual System::Collections::Generic::Dictionary< int, array< int16_t >^> ^ GetGroupChannelDatal16 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< uint32_t >^> ^ GetGroupChannelDataUl32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

• virtual System::Collections::Generic::Dictionary< int, array< int32_t $>^{\land}>$ $^{\land}$ GetGroupChannelDatal32 (DacqGroupChannelEnumNet group, int frames, [System::Runtime::InteropServices::Out]int % frames_ret)

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

- void SetupGroupDacqQueue (int queuesize, int threshold)
- void SetupGroupDacqQueue (int queuesize, int threshold, unsigned int virtualDevice)
- CHWInfo ^ HWInfo ()

Static Public Attributes

- static const int Error_Callback_Queue_Full = 0x100
- static const int Error_Callback_Aquisition_Stopped = 0x200
- static const int Error Callback Packet Error = 1
- static const int Error_Callback_RingQueue_Full = 3
- static const int Error Callback Frames Lost = 4
- static const int Error_Callback_Data_lost = 5

Properties

virtual int Samplerate [get, set]
 The sampling frequency of the device in Hz.

Events

- OnChannelData^ ChannelDataEvent [add, remove, raise]
- OnError[∧] ErrorEvent [add, remove, raise]

Additional Inherited Members

11.52.1 Detailed Description

Base class for data acquisition devices.

11.52.2 Constructor & Destructor Documentation

```
11.52.2.1 CMcsUsbDacqNet() CMcsUsbDacqNet ()
```

```
11.52.2.2 ~ CMcsUsbDacqNet() ~ CMcsUsbDacqNet ( )
```

11.52.3 Member Function Documentation

```
11.52.3.1 AddSelectedChannelsQueue() [1/4] virtual int AddSelectedChannelsQueue (
```

```
int nByteOffset,
int nChannelOffset,
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize ) [virtual]
```

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read← FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

Parameters

Parameters

queuesize	Size of sample frames the FIFO can hold.
-----------	--

Parameters

threshold	Number of sample frames the FIFO must ac	quire before the callback function is called.
-----------	--	---

Parameters

samplesize size of the datawords, either 16 or 32bit
--

Returns

The handle to the Queue.

$\textbf{11.52.3.2} \quad \textbf{AddSelectedChannelsQueue() [2/4]} \quad \text{virtual int AddSelectedChannelsQueue (} \\$

```
int nByteOffset, int nChannelOffset, array< bool >^{\wedge} selectedChannels,
```

```
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize ) [virtual]
```

11.52.3.3 AddSelectedChannelsQueue() [3/4] virtual int AddSelectedChannelsQueue (

```
int nByteOffset,
int nChannelOffset,
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize ) [virtual]
```

Adds a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_Read FramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

nByteOffset Nur	nber of bytes to start with.
-----------------	------------------------------

Parameters

1	nChannelOffset	Number of channel to start with (counted in samplesize bytes).	
---	----------------	--	--

Parameters

nChannels	Number of channels to be collected in the FIFO.
Honamicis	i Number of Chambers to be conected in the fin O.

Parameters

ueuesize Size of sample frames the FIFO can	nold.
---	-------

threshold Number of sample frames the FIFO must acquire before the callback function is cal	threshold
---	-----------

samplesize	size of the datawords, either 16 or 32bit.
------------	--

Returns

The handle to the Queue.

```
11.52.3.4 AddSelectedChannelsQueue() [4/4] virtual int AddSelectedChannelsQueue (
```

```
int nByteOffset,
int nChannelOffset,
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize) [virtual]
```

```
11.52.3.5 ChannelBlock_AvailFrames() [1/2] virtual uint32_t ChannelBlock_AvailFrames ( int handle ) [virtual]
```

Get the number of sample frames already available in the FIFO.

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.

Returns

Number of sample frames available in the FIFO.

```
11.52.3.6 ChannelBlock_AvailFrames() [2/2] virtual uint32_t ChannelBlock_AvailFrames (
    int handle,
    int queue ) [virtual]
```


[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Parameters

frames	Number of sample frames to read.
--------	----------------------------------

Parameters

```
frames_ret Number of sample frames which were read, might be smaller than frames.
```

Returns

Array of int16_t frame arrays.

```
11.52.3.8 ChannelBlock_ReadAsFrameArrayl16() [2/2] virtual array<array<int16_t>^> ^ Channel←
Block_ReadAsFrameArrayI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

```
handle Handle of the FIFO queue. Zero when the SetSelectedData call was used.
```

queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

	frames_ret	Number of sample frames which were read, might be smaller than frames.
--	------------	--

Returns

Array of int16_t frame arrays.

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Parameters

frames	Number of sample frames to read.
--------	----------------------------------

Returns

Array of int32_t frame arrays.

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

Returns

Array of int32_t frame arrays.

Read data from a FIFO queue as array of uint16_t data frame arrays

handle Handle of the FIFO queue. Zero when the SetSelectedData call was used
--

frames Number of sample frames to read.

Parameters

Returns

Array of uint16_t frame arrays.

$\textbf{11.52.3.12} \quad \textbf{ChannelBlock_ReadAsFrameArrayUl16() [2/2]} \quad \text{virtual array} < \text{array} < \text{uint16_t} > ^ \\ > \ ^ \\ \\$

```
ChannelBlock_ReadAsFrameArrayUI16 (
    int handle,
    int queue,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedData call was used.
--------	--

Parameters

queue	Number of the sub queue.
frames	Number of sample frames to read.

were read, might be smaller than frames.
--

Returns

Array of uint16_t frame arrays.

Read data from a FIFO queue as array of uint16_t data frame arrays

Parameters

Parameters

Number of sample frames to r	ead.
------------------------------	------

Parameters

```
frames_ret Number of sample frames which were read, might be smaller than frames.
```

Returns

Array of uint32_t frame arrays.


```
int nandle,
int queue,
int frames,
[System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue as array of uint16_t data frame arrays

handle	Handle of the FIFO queue. Zero when the SetSelectedData call was used.

Parameters

queue	Number of the sub queue.
frames	Number of sample frames to read.

Parameters

Returns

Array of uint32_t frame arrays.

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

, ,,	11 1 1 5 5 1 1 0 0 1 1 0 1 1 1	
nangie i	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.	
	That are of the first queen and the first are detected and the detect	

frames Number of sample frames to read.	Number of sample frames to read.	
	frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int16_t arrays and hardware channel as key.

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

handle Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call v	vas used.
--	-----------

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int32_t arrays and hardware channel as key.

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

```
handle Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
```

frames Number of sample frames to read.	
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint16_t arrays and hardware channel as key.

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

handle	Handle of the FIFO queue. Zero when the SetSelectedChannelsQueue call was used.
--------	---

Parameters

fr	rames	Number of sample frames to read.
fr	rames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint32_t arrays and hardware channel as key.

```
11.52.3.19 ChannelBlock_ReadFramesI16() [1/2] virtual void ChannelBlock_ReadFramesI16 (
    int handle,
    array< int16_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in int16_t data format

handle	Handle of the FIFO queue.	Either zero when the SetSelectedData call was used or the channel number.
--------	---------------------------	---

Parameters

buffer	Buffer to put the data from the device in.
frames_pos	Position in buffer where to put the data.
frames	Number of sample frames to read.

Parameters

Read data from a FIFO queue in int16_t data format

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the cha	nnel number.
--	--------------

Parameters

frames Number of sample frames to read
--

ſ		
١	trames ret	Number of sample frames which were read, might be smaller than frames.
۱		rtambor or bampio marrico trinori troto roda, migne bo omanor triar marrico.

```
11.52.3.21 ChannelBlock_ReadFramesI32() [1/2] virtual void ChannelBlock_ReadFramesI32 (
    int handle,
    array< int32_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

queue. Either zero when the SetSelectedData call was used or the channel number.	handle
--	--------

Parameters

buffer	Buffer to put the data from the device in.
frames_pos	Position in buffer where to put the data.
frames	Number of sample frames to read.

Parameters

ret Number of sample frames which were read, might be smaller than frames.
--

```
11.52.3.22 ChannelBlock_ReadFramesl32() [2/2] virtual array<int32_t> ^ ChannelBlock_Read←
FramesI32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

handle Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.

<i>trames</i> Number of sample frames to read.	frames	Number of sample frames to read.
--	--------	----------------------------------

Parameters

	frames_ret	Number of sample frames which were read, might be smaller than frames.
--	------------	--

```
11.52.3.23 ChannelBlock_ReadFramesUl16() [1/2] virtual void ChannelBlock_ReadFramesUl16 (
    int handle,
    array< uint16_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint16_t data format

Parameters

handle	Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.	
	number.	
buffer	Buffer to put the data from the device in.	
frames_pos	Position in buffer where to put the data.	
frames	Number of sample frames to read.	

Parameters

```
frames_ret Number of sample frames which were read, might be smaller than frames.
```

```
11.52.3.24 ChannelBlock_ReadFramesUl16() [2/2] virtual array<uint16_t> ^ ChannelBlock_Read←
FramesUI16 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint16_t data format

handle Handle of the FIFO gueue. Either zero when the SetSelectedData call was used or the channel number

Parameters

Parameters

	frames_ret	Number of sample frames which were read, might be smaller than frames.
--	------------	--

Returns

Array of data from the device.

```
11.52.3.25 ChannelBlock_ReadFramesUl32() [1/2] virtual void ChannelBlock_ReadFramesUI32 (
    int handle,
    array< uint32_t >^ buffer,
    int frames_pos,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

buffer	Buffer to put the data from the device in.
frames_pos	Position in buffer where to put the data.
frames	Number of sample frames to read.

frames_ret | Number of sample frames which were read, might be smaller than frames.

```
11.52.3.26 ChannelBlock_ReadFramesUl32() [2/2] virtual array<uint32_t> ^ ChannelBlock_Read←
FramesUl32 (
    int handle,
    int frames,
    [System::Runtime::InteropServices::Out] int % frames_ret ) [virtual]
```

Read data from a FIFO queue in uint32_t data format

Parameters

handle | Handle of the FIFO queue. Either zero when the SetSelectedData call was used or the channel number.

Parameters

frames | Number of sample frames to read.

Parameters

frames_ret | Number of sample frames which were read, might be smaller than frames.

11.52.3.27 ClearBuffers() virtual void ClearBuffers () [virtual]

```
11.52.3.29 GetAdapterType() virtual AdapterTypeEnumNet GetAdapterType ( ) [virtual]
```

Gets the adapter which is connected to the MEA2100 device.

Returns

AdapterTypeEnumNet which enumerates the possible adapters.

```
11.52.3.30 GetAdcDataFormat() virtual uint32_t GetAdcDataFormat ( uint32_t virtualDevice ) [virtual]
```

Gets the ADC data format, 16 means 16 bits, 24 means 24 bits, 32 means 32 bits.

Returns

The data format in bits.

11.52.3.33 GetChannelDataFillSize() int GetChannelDataFillSize ()

```
11.52.3.35 GetChannelsInBlock() virtual int32_t GetChannelsInBlock ( unsigned int virtualDevice ) [virtual]
```

Get the number of 16 bit datawords which will be collected per sample frame, use after the device is configured.

Returns

Number of 16 bit datawords per sample frame.

```
11.52.3.37 GetDataMode() virtual DataModeEnumNet GetDataMode ( unsigned int virtualDevice ) [virtual]
```

Gets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

virtualDevice Virtual device to use

Returns

DataModeEnumNet which enumerates the possible data modes.

Gets the function/source of an digital output bit.

This overload is for the CMOSMEA5000 device.

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This is the templated generic implementation.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the MEA2100 device.

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the Multiwell device.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Gets the function/source of an digital output bit.

This overload is for the SCU device.

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

```
11.52.3.44 GetDigitalSource() [7/7] void GetDigitalSource (

DigitalTargetEnumNet digitaltarget,
```

```
int32_t NrChannel,
[System::Runtime::InteropServices::Out] W2100DigitalSourceEnumNet% source,
[System::Runtime::InteropServices::Out] int% bitnumber_offset )
```

Gets the function/source of an digital output bit.

This overload is for the W2100 device.

Parameters

digitaltarget	The digital target to query.
NrChannel	The channel/bit of target to query.
source	The source/function assignd to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

11.52.3.45 GetErrorMessage() uint32_t GetErrorMessage (

DacqGroupChannelEnumNet group,

Read data from a FIFO queue in int16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

[System::Runtime::InteropServices::Out] int % frames_ret) [virtual]

Parameters

grou	ıp	Group selector supported by the device.
------	----	---

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int16_t arrays and hardware channel as key.

Read data from a FIFO queue in int32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

group	Group selector supported by the device.
-------	---

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of int32_t arrays and hardware channel as key.

Read data from a FIFO queue in uint16_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

aroun	Group selector supported by the device.

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint16_t arrays and hardware channel as key.

Read data from a FIFO queue in uint32_t data format, that contains subqueues, each populates an entry in the dictionary by hardware channel number

Parameters

```
group Group selector supported by the device.
```

Parameters

frames	Number of sample frames to read.
frames_ret	Number of sample frames which were read, might be smaller than frames.

Returns

Dictonary of uint32_t arrays and hardware channel as key.

Gets the maximal sampling frequency of the device.

Returns

Sampling frequency in Hz.

```
11.52.3.54 GetMeaLayout() virtual MeaLayoutEnumNet GetMeaLayout ( ) [virtual]
```

Gets the MEA layout which is connected to the MEA2100 device.

Returns

MeaLayoutEnumNet which enumerates the MEA types.

```
11.52.3.55 GetMinSamplingFrequencyStepsize() virtual uint32_t GetMinSamplingFrequencyStepsize ( ) [virtual]
```

Gets the minimal sampling frequency step size increment value of the device.

Returns

Sampling frequency step size in Hz.

Get the real number of data bits.

This value may be different from the value returned by GetDataFormat, e.g. in MC_Card the data are shifted 2 bits so the real number is 14 while the data format is 16 bits

```
11.52.3.57 GetPoti() virtual uint32_t GetPoti (
             uint32_t channel,
              [System::Runtime::InteropServices::Out] \ uint32\_t\% \ value \ ) \quad [virtual]
11.52.3.58 GetResolutionPerDigit() virtual uint32_t GetResolutionPerDigit (
             uint32_t virtualDevice,
             DacqGroupChannelEnumNet group,
             [System::Runtime::InteropServices::Out] int% res,
              [System::Runtime::InteropServices::Out] int% resUnit ) [virtual]
11.52.3.59 GetSamplerate() virtual int32_t GetSamplerate (
             unsigned int virtualDevice ) [virtual]
Gets the sampling frequency of the device.
Returns
     Sampling frequency in Hz.
11.52.3.60 GetVoltageRangeIndex() virtual uint32_t GetVoltageRangeIndex (
             unsigned int virtualDevice ) [virtual]
11.52.3.61 GetVoltageRangeInMicroVolt() virtual int32_t GetVoltageRangeInMicroVolt (
             unsigned int virtualDevice ) [virtual]
Gets the currently selected voltage range on devices which support multiple voltage ranges.
Returns
     The Voltage Range in uV.
11.52.3.62 GetVoltageRangeInMilliVolt() virtual int32_t GetVoltageRangeInMilliVolt () [virtual]
```

Returns

The rounded Voltage Range in mV.

Gets the currently selected voltage range on devices which support multiple voltage ranges.

```
11.52.3.63 HWInfo() CHWInfo ^ HWInfo ( )

11.52.3.64 SendStartDacq() [1/2] virtual void SendStartDacq ( ) [virtual]

Start sampling.

11.52.3.65 SendStartDacq() [2/2] virtual void SendStartDacq ( int VirtualDacqMap ) [virtual]

Start sampling.

Parameters
```

Start sampling together with the STG.

Parameters

trigger_map
VirtualDacqMap

VirtualDacqMap

11.52.3.67 SendStopDacq() [1/2] virtual void SendStopDacq () [virtual]

Stop sampling.

```
11.52.3.68 SendStopDacq() [2/2] virtual void SendStopDacq ( int VirtualDacqMap ) [virtual]
```

Stop sampling.

D					
Pa	ra	m	ല	ſΑ	rς

VirtualDacqMap

Stop sampling together with the STG.

Parameters

trigger_map

Stop sampling together with the STG and options.

Parameters

trigger_map

Parameters

options

Parameters

VirtualDacqMap

```
11.52.3.71 SetDataMode() virtual void SetDataMode (

DataModeEnumNet dataMode,

unsigned int virtualDevice ) [virtual]
```

Sets the data mode, can be 16, 24 or 32bit, all signed or unsigned on the MEA2100 device.

Parameters

dataMode	DataModeEnumNet enumerates the possible data modes.
virtualDevice	Virtual device to use.

Sets the function/source of an digital output bit.

This overload is for the CMOSMEA5000 device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This is the templated generic implementation.

	digitaltarget	The digital target to change.
ſ	NrChannel	The channel/bit of target to change.
Ī	source	The source/function to assign to the digital target.
Ī	bitnumber offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the MEA2100-256 device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the MEA2100 device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the Multiwell device.

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the SCU device.

Parameters

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

Sets the function/source of an digital output bit.

This overload is for the W2100 device.

digitaltarget	The digital target to change.
NrChannel	The channel/bit of target to change.
source	The source/function to assign to the digital target.
bitnumber_offset	An offset / bit number with the source/function.

```
11.52.3.80 SetSamplerate() virtual void SetSamplerate ( int32_t rate,
```

```
unsigned int oversample,
unsigned int virtualDevice ) [virtual]
```

Sets the sampling frequency of the device.

```
rate Sampling frequency in Hz.
```


Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be devided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

selectedChannels	List of channels to be collected in the FIFO.
------------------	---

Parameters

ple frames the FIFO can hold.	queuesize Size of sample
-------------------------------	----------------------------

Parameters

threshold	Number of sample frames the FIFO must acquire before the callback function is called.
-----------	---

Parameters

samplesize	size of the datawords, either 16 or 32bit.
------------	--

ChannelsInBlock	value obtained from GetChannelsInBlock.
-----------------	---

int ChannelsInBlock) [virtual]

```
11.52.3.83 SetSelectedChannels() [3/4] virtual void SetSelectedChannels (
    int nChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    int ChannelsInBlock ) [virtual]
```

Create a FIFO queue per channel. Each channel will have its own FIFO and Callback function.

When using a 32bit sample size, the number obtained from GetChannelsInBlock must be devided by 2 to be used here, since GetChannelsInBlock returns the number of 16 bit datapoints per sample frame, while this functions uses the number of sample frames in its own data format.

Parameters

Parameters

Parameters

mber of samples frames the FIFO must acquire before the callback function is called	
---	--

samplesize	size of the datawords, either 16 or 32bit.
ChannelsInBlock	value obtained from GetChannelsInBlock.

11.52.3.84 SetSelectedChannels() [4/4] virtual void SetSelectedChannels (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

11.52.3.85 SetSelectedChannelsQueue() [1/4] virtual void SetSelectedChannelsQueue (

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock) [virtual]
```

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

selectedChannels	List of channels to be collected in the FIFO.
SCICCICUCITALITICIS	LISE OF CHAILIEIS TO DE COHECTEU III THE FIFO.

Parameters

queuesize	Size of sample frames the FIFO can hold.
queuesize	Size of sample frames the Fit O can hold

Parameters

threshold	Number of sample frames the FIFO must acquire before the callba	ack function is called.
-----------	---	-------------------------

campleciza	size of the datawords, either 16 or 32bit.
Sallibicsize	Size of the datawords, ettiler to or szbit.

```
11.52.3.87 SetSelectedChannelsQueue() [3/4] virtual void SetSelectedChannelsQueue (
    int nChannels,
    int queuesize,
```

int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock) [virtual]

 $\verb|int| ChannelsInBlock|) | [virtual]$

Create a common FIFO queue for all channels. Data in callback will be a list per channel. Use ChannelBlock_← ReadFramesDict... with handle = 0 to read the data.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

nChannels	Number of channels to be collected in the FIFO.
-----------	---

Parameters

queuesize	Size of sample frames the FIFO can hold.
-----------	--

threshold Number of sample frames the FIFO must acquire before the callback function is called.

samplesize	size of the datawords, either 16 or 32bit.
------------	--

Parameters

11.52.3.88 SetSelectedChannelsQueue() [4/4] virtual void SetSelectedChannelsQueue (

```
int nChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
SampleDstSizeNet sampleDstSize,
int ChannelsInBlock ) [virtual]
```

11.52.3.89 SetSelectedData() [1/4] virtual void SetSelectedData (

```
array< bool >^ selectedChannels,
int queuesize,
int threshold,
SampleSizeNet samplesize,
int ChannelsInBlock ) [virtual]
```

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

selectedChannels	List of channels to be collected in the FIFO.
------------------	---

aueuesize	Size of sample frames the FIFO can hold.
446463126	, Olze of sample marries the rin o can hold.

Parameters

samplesize	size of the datawords, either 16 or 32bit.
ChannelsInBlock	value obtained from GetChannelsInBlock.

```
11.52.3.91 SetSelectedData() [3/4] virtual void SetSelectedData (
    int nChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    int ChannelsInBlock ) [virtual]
```

Create a common FIFO queue for all channels. Use handle = 0 in the ChannelBlock_ReadFrames... functions.

When using 32 bit data format, ChannelsInBlock is still the number of 16 bit channels per frame, as obtained from GetChannelsInBlock, while nChannels is the number of 32 bit channels to be read from the device. So when all channels from a device are read in 32 bit data format nChannels = ChannelsInBlock/2

Parameters

	nChannels	Number of channels to be collected in the FIFO.
--	-----------	---

queuesize Size of sample	e frames the FIFO can hold.
--------------------------	-----------------------------

thre	shold	Number of sample frames the FIFO must acquire before the callback function is called.
------	-------	---

Parameters

samplesize	size of the datawords, either 16 or 32bit.	
ChannelsInBlock	value obtained from GetChannelsInBlock.	

```
11.52.3.92 SetSelectedData() [4/4] virtual void SetSelectedData (
    int nChannels,
    int queuesize,
    int threshold,
    SampleSizeNet samplesize,
    SampleDstSizeNet sampleDstSize,
    int ChannelsInBlock ) [virtual]

11.52.3.93 SetupGroupDacqQueue() [1/2] void SetupGroupDacqQueue (
    int queuesize,
    int threshold )

11.52.3.94 SetupGroupDacqQueue() [2/2] void SetupGroupDacqQueue (
    int queuesize,
    int threshold,
    unsigned int virtualDevice )
```

Sets the voltage range on devices which support multiple voltage ranges.

unsigned int virtualDevice) [virtual]

int32_t voltageRangeIndex,

Parameters

 $\textbf{11.52.3.95} \quad \textbf{SetVoltageRangeByIndex()} \quad \texttt{virtual void SetVoltageRangeByIndex} \ \ \textbf{(}$

Sets the voltage range on devices which support multiple voltage ranges.

Parameters

voltageRange Voltage Range to use in μV

This replaces SetVoltageRange, where the value of the range was in mV!

```
11.52.3.97 StartDacq() [1/4] virtual void StartDacq ( ) [virtual]
```

Start the data acquisition thread and sampling.

```
11.52.3.98 StartDacq() [2/4] virtual void StartDacq ( int32_t timeout ) [virtual]
```

Start the data acquisition thread and sampling.

Parameters

```
timeout Timeout in ms.
```

Start the data acquisition thread and sampling.

timeout	Timeout in ms.
---------	----------------

numSubmittedUsbBuffers	Number of USB Buffers that are simultaniously submitted.	

Parameters

Parameters

packetsInUrb	Packets in each URB.
--------------	----------------------

Start the data acquisition thread and sampling.

Parameters

numSubmittedUsbBuffers	Number of USB Buffers that are simultaniously submitted.

Parameters

timeout	Timeout in ms.
unicoul	i iiiiicout iii iiio.

Parameters

11.52.3.101 StartLoop() [1/4] virtual void StartLoop () [virtual]

Start the data acquisition thread.

```
11.52.3.102 StartLoop() [2/4] virtual void StartLoop ( int32_t timeout ) [virtual]
```

Start the data acquisition thread.

Parameters

```
timeout Timeout in ms.
```

Start the data acquisition thread.

Parameters

timeout	Timeout in ms.
umeoui	i illieout ili ilis.

numSubmittedUsbBuffers	Number of USB Buffers that are simultaniously submitted.

numUsbBuffers	Number of USB Buffers to use.
---------------	-------------------------------

Parameters

packetsInUrb | Packets in each URB.

Start the data acquisition thread.

Parameters

numSubmittedUsbBuffers	Number of USB Buffers that are simultaniously submitted.
------------------------	--

Parameters

timeout	Timeout in ms.
---------	----------------

Parameters

numUsbButters Nui	nber of USB Buffers to use.
---------------------	-----------------------------

packetsInUrb Packets in each URB.

11.52.3.105 StopDacq() [1/2] virtual void StopDacq () [virtual]

Stop the data acquisition thread and sampling.

11.52.3.106 StopDacq() [2/2] virtual void StopDacq (uint32_t virtualDevice) [virtual]

Stop the data acquisition thread and sampling.

Parameters

virtualDevice	Virtual Device to start.
---------------	--------------------------

11.52.3.107 StopLoop() virtual void StopLoop () [virtual]

11.52.4 Member Data Documentation

11.52.4.1 Error_Callback_Aquisition_Stopped const int Error_Callback_Aquisition_Stopped = 0x200 [static]

11.52.4.2 Error_Callback_Data_lost const int Error_Callback_Data_lost = 5 [static]

11.52.4.3 Error_Callback_Frames_Lost const int Error_Callback_Frames_Lost = 4 [static]

11.52.4.4 Error_Callback_Packet_Error const int Error_Callback_Packet_Error = 1 [static]

11.52.4.5 Error_Callback_Queue_Full const int Error_Callback_Queue_Full = 0x100 [static]

11.52.4.6 Error_Callback_RingQueue_Full const int Error_Callback_RingQueue_Full = 3 [static]

11.52.5 Property Documentation

11.52.5.1 Samplerate virtual int Samplerate [get], [set]

The sampling frequency of the device in Hz.


11.52.6 Event Documentation

11.52.6.1 ChannelDataEvent OnChannelData^ ChannelDataEvent [add], [remove], [raise]

11.52.6.2 ErrorEvent OnError^ ErrorEvent [add], [remove], [raise]

11.53 CMcsUsbDeviceStatePushFunctionNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushFunctionNet:

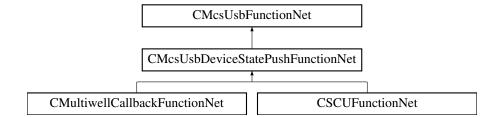

Public Member Functions

void TriggerStatus ()

Protected Member Functions

CMcsUsbDeviceStatePushFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] p
 Device)

Events

• OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add, remove, raise]

Additional Inherited Members

11.53.1 Constructor & Destructor Documentation

```
11.53.1.1 CMcsUsbDeviceStatePushFunctionNet() CMcsUsbDeviceStatePushFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pDevice ) [protected]
```

11.53.2 Member Function Documentation


```
11.53.2.1 TriggerStatus() void TriggerStatus ( )
```

11.53.3 Event Documentation

11.53.3.1 McsUsbDeviceStateEvent OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]

11.54 CMcsUsbDeviceStatePushNet Class Reference

Inheritance diagram for CMcsUsbDeviceStatePushNet:

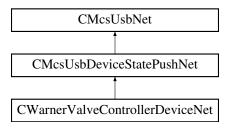

Public Member Functions

• void TriggerStatus ()

Protected Member Functions

CMcsUsbDeviceStatePushNet (CMcsUsbPointerContainer[^] pDevice)

Events

• OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add, remove, raise]

Additional Inherited Members

11.54.1 Constructor & Destructor Documentation

```
11.54.1.1 CMcsUsbDeviceStatePushNet() CMcsUsbDeviceStatePushNet (
CMcsUsbPointerContainer^ pDevice) [protected]
```

11.54.2 Member Function Documentation


```
11.54.2.1 TriggerStatus() void TriggerStatus ( )
```

11.54.3 Event Documentation

```
11.54.3.1 McsUsbDeviceStateEvent OnMcsUsbDeviceState^ McsUsbDeviceStateEvent [add], [remove], [raise]
```

11.55 CMcsUsbFactoryNet Class Reference

Inheritance diagram for CMcsUsbFactoryNet:

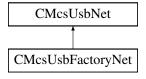

Classes

· class TargetSelectionInfoNet

Public Member Functions

- CMcsUsbFactoryNet ()
- ∼CMcsUsbFactoryNet ()
- unsigned int GetNumDestinations ()
- String \(^\) GetDestinationName (unsigned int index)
- String \(^\) GetDestinationName (CFirmwareDestinationNet dest)
- void SetDestinationSerialNumber (CFirmwareDestinationNet dest, String[^] serialnumber)
- String ^ GetDestinationSerialNumber (CFirmwareDestinationNet dest)
- CFirmwareDestinationNet GetDestination (unsigned int index)
- CFirmwareDestinationNet GetDestination (String[^] Key)
- unsigned int GetDestinationTargetAddress (CFirmwareDestinationNet destination)

Gets the target base address for the destination.

- TargetSelectionInfoNet ^ GetDestinationTargetSelectionInfo (CFirmwareDestinationNet destination)
- uint32 t ChangeSerialNumber (String^ serial)
- bool LoadUserFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry)

Send the DSP Firmware to the MEA21 device.

- bool LoadUserFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, uint32_t LockMask)
- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, bool SkipWait)

Flashes a firmware file to the device.

- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, bool SkipWait, unsigned int LockMask)
- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest)

Flashes a firmware file to the device.

 bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest, bool SkipWait)

Flashes a firmware file to the device.

- bool UpdateFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet dest, bool SkipWait, uint32 t LockMask)
- bool CompareFirmware (String[^] FirmwareFile, CMcsUsbListEntryNet[^] listEntry, CFirmwareDestinationNet Dest, OnUpdateFirmwareStatusChange[^] deleg, OnUpdateFirmwareProgress[^] progress, String[^] MessagePrefix, unsigned int LockMask, [System::Runtime::InteropServices::Out] String[^]% ErrorText, [System::Runtime::InteropServices::Out] String[^]% Protokoll)
- uint32 t Coldstart (CFirmwareDestinationNet dest)
- int32 t GetXilinxFlashOffset (CFirmwareDestinationNet dest)
- uint32_t GetXilinxFlashReadCommand (CFirmwareDestinationNet dest)
- array< uint8_t > ^ DownloadFirmware (CFirmwareDestinationNet Dest, uint32_t Address, uint32_t length)
- array< unsigned char > ^ ReadBlockFromFlash (unsigned int FPGA, unsigned int Address)
- void ReadBlockFromFlash (unsigned int FPGA, unsigned int Address, array< unsigned char >^ buffer, int position)
- array< unsigned char > ^ ReadBlockFromIFBGlobalEEprom (unsigned int Address)
- array< unsigned char > ^ ReadBlockFromNVMEM (unsigned int FPGA, unsigned int Offset, unsigned int Address)

Static Public Member Functions

- static String ^ GetDestinationDisplayLabel (String^ RawLabel, CFirmwareDestinationNet dest)
- static String [^] FindFirmwareVersionMagicInBuffer (array< unsigned char > ^ buffer, int length, [System::
 — Runtime::InteropServices::Out]int% position)
- static bool GetFirmwareVersionFromFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)

Retrives version info from a Firmware update file.

- static bool GetFirmwareVersionFromFile (String^ FirmwareFile, [System::Runtime::InteropServices::Out] uint32 t% Version, [System::Runtime::InteropServices::Out] uint32 t% Position)
- static bool GetFirmwareVersionFromHexFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] uint32_t% Version)
- static uint32_t GetChecksumFromFX3Image (String[^] FirmwareFile)
- static uint32_t GetUSBDeviceIDFromFX3Image (String^{\(\)} FirmwareFile)
- static bool GetUsercodeFromBitFile (String[^] FirmwareFile, [System::Runtime::InteropServices::Out] unsigned int% Usercode)

Static Public Attributes

- static const uint32_t FX3MCSDataAddress = 0x40037E00
- static const uint32 t FX3MCSDataDeviceIdOffset = 0x4
- static const uint32 t FX3MCSDataVersionOffset = 0x8
- static const uint32_t FX3MCSDatalFB2ImageOffset = 0xC
- static const uint32_t FX3MCSDataIFB1ImageOffset = 0x2C

Additional Inherited Members

11.55.1 Constructor & Destructor Documentation

```
11.55.1.1 CMcsUsbFactoryNet() CMcsUsbFactoryNet ( )
```

```
11.55.1.2 ~CMcsUsbFactoryNet() ~CMcsUsbFactoryNet ()
```

11.55.2 Member Function Documentation

```
11.55.2.1 ChangeSerialNumber() uint32_t ChangeSerialNumber ( String^ serial)
```

```
11.55.2.2 Coldstart() uint32_t Coldstart (
             CFirmwareDestinationNet dest )
11.55.2.3 CompareFirmware() bool CompareFirmware (
             String^ FirmwareFile,
             CMcsUsbListEntryNet^ listEntry,
             CFirmwareDestinationNet Dest,
             OnUpdateFirmwareStatusChange^ deleg,
             OnUpdateFirmwareProgress^ progress,
             String Message Prefix,
             unsigned int LockMask,
             [System::Runtime::InteropServices::Out] String^{\%} ErrorText,
             [System::Runtime::InteropServices::Out] \ String^{\ } \ \textit{Protokoll} \ )
11.55.2.4 DownloadFirmware() array<uint8_t> ^ DownloadFirmware (
             CFirmwareDestinationNet Dest,
             uint32_t Address,
             uint32_t length )
11.55.2.5 FindFirmwareVersionMagicInBuffer() static String ^ FindFirmwareVersionMagicInBuffer (
             array< unsigned char >^{\wedge} buffer,
             int length,
             [System::Runtime::InteropServices::Out] int% position ) [static]
11.55.2.6 GetChecksumFromFX3Image() static uint32_t GetChecksumFromFX3Image (
             String<sup>∧</sup> FirmwareFile ) [static]
11.55.2.7 GetDestination() [1/2] CFirmwareDestinationNet GetDestination (
             String^ Key )
11.55.2.8 GetDestination() [2/2] CFirmwareDestinationNet GetDestination (
             unsigned int index)
11.55.2.9 GetDestinationDisplayLabel() static String ^ GetDestinationDisplayLabel (
             String^ RawLabel,
             CFirmwareDestinationNet dest ) [static]
```

```
11.55.2.10 GetDestinationName() [1/2] String ^ GetDestinationName (
              CFirmwareDestinationNet dest )
11.55.2.11 GetDestinationName() [2/2] String ^ GetDestinationName (
              unsigned int index )
11.55.2.12 GetDestinationSerialNumber() String ^ GetDestinationSerialNumber (
              CFirmwareDestinationNet dest )
11.55.2.13 GetDestinationTargetAddress() unsigned int GetDestinationTargetAddress (
              CFirmwareDestinationNet destination )
Gets the target base address for the destination.
Parameters
 destination
              The destination to be queried.
Returns
     The base address as a 32 bit number, only the lower 16 bit represent the address.
11.55.2.14 GetDestinationTargetSelectionInfo() TargetSelectionInfoNet ^ GetDestinationTarget←
SelectionInfo (
              CFirmwareDestinationNet destination )
\textbf{11.55.2.15} \quad \textbf{GetFirmwareVersionFromFile() [1/2]} \quad \texttt{static bool GetFirmwareVersionFromFile ()} \\
              String^ FirmwareFile,
              [System::Runtime::InteropServices::Out] uint32_t% Version ) [static]
Retrives version info from a Firmware update file.
11.55.2.16 GetFirmwareVersionFromFile() [2/2] static bool GetFirmwareVersionFromFile (
              String Firmware File,
              [System::Runtime::InteropServices::Out] uint32_t% Version,
```

[System::Runtime::InteropServices::Out] uint32_t% Position) [static]

```
11.55.2.17 GetFirmwareVersionFromHexFile() static bool GetFirmwareVersionFromHexFile (
             String^ FirmwareFile,
             [System::Runtime::InteropServices::Out] uint32_t% Version ) [static]
11.55.2.18 GetNumDestinations() unsigned int GetNumDestinations ( )
11.55.2.19 GetUSBDeviceIDFromFX3Image() static uint32_t GetUSBDeviceIDFromFX3Image (
             String FirmwareFile ) [static]
11.55.2.20 GetUsercodeFromBitFile() static bool GetUsercodeFromBitFile (
             String^ FirmwareFile,
             [System::Runtime::InteropServices::Out] unsigned int% Usercode ) [static]
11.55.2.21 GetUsercodeFromFlash() bool GetUsercodeFromFlash (
             unsigned int FPGA,
             unsigned int Address,
             [System::Runtime::InteropServices::Out] unsigned int% Usercode )
11.55.2.22 GetXilinxFlashOffset() int32_t GetXilinxFlashOffset (
             CFirmwareDestinationNet dest )
11.55.2.23 GetXilinxFlashReadCommand() uint32_t GetXilinxFlashReadCommand (
             CFirmwareDestinationNet dest )
11.55.2.24 LoadUserFirmware() [1/2] bool LoadUserFirmware (
             String^ FirmwareFile,
             CMcsUsbListEntryNet^ listEntry )
Send the DSP Firmware to the MEA21 device.
Parameters
 FirmwareFile
              Filename of the DSP Firmware (*.bin) file.
```

listEntry Device to use for the connection. See CMcsUsbListNet.

```
11.55.2.25 LoadUserFirmware() [2/2] bool LoadUserFirmware (
               String^ FirmwareFile,
               CMcsUsbListEntryNet<sup>^</sup> listEntry,
               uint32_t LockMask )
\textbf{11.55.2.26} \quad \textbf{ReadBlockFromFlash() [1/2]} \quad \text{array} < \text{unsigned char} > \ ^ \land \quad \textbf{ReadBlockFromFlash ()}
               unsigned int FPGA,
               unsigned int Address )
11.55.2.27 ReadBlockFromFlash() [2/2] void ReadBlockFromFlash (
               unsigned int FPGA,
               unsigned int Address,
               array< unsigned char >^{\wedge} buffer,
               int position )
11.55.2.28 ReadBlockFromIFBGlobalEEprom() array<unsigned char> ^ ReadBlockFromIFBGlobal←
EEprom (
               unsigned int Address )
\textbf{11.55.2.29} \quad \textbf{ReadBlockFromNVMEM()} \quad \texttt{array} < \texttt{unsigned char} > \ ^ \land \ \texttt{ReadBlockFromNVMEM} \ (
               unsigned int FPGA,
               unsigned int Offset,
               unsigned int Address )
11.55.2.30 SetDestinationSerialNumber() void SetDestinationSerialNumber (
               CFirmwareDestinationNet dest,
               String^{\wedge} serialnumber)
11.55.2.31 UpdateFirmware() [1/5] bool UpdateFirmware (
               String^{\wedge} FirmwareFile,
               CMcsUsbListEntryNet^ listEntry,
               CFirmwareDestinationNet dest )
```

Flashes a firmware file to the device.

FirmwareFile Filename of the Firmware file.	∍. │
---	------

Parameters

```
listEntry Device to use for the connection.
```

Flashes a firmware file to the device.

Parameters

FirmwareFile File	name of the Firmware file.
-------------------	----------------------------

```
listEntry Device to use for the connection.
```

```
OnUpdateFirmwareStatusChange^ deleg,
OnUpdateFirmwareProgress^ progress,
bool SkipWait)
```

Flashes a firmware file to the device.

Parameters

```
FirmwareFile | Filename of the Firmware file.
```

```
11.55.2.35 UpdateFirmware() [5/5] bool UpdateFirmware (
```

```
String FirmwareFile,

CMcsUsbListEntryNet listEntry,

CFirmwareDestinationNet Dest,

OnUpdateFirmwareStatusChange deleg,

OnUpdateFirmwareProgress progress,

bool SkipWait,

unsigned int LockMask)
```


11.55.3 Member Data Documentation

```
11.55.3.1 FX3MCSDataAddress const uint32_t FX3MCSDataAddress = 0x40037E00 [static]
```

- 11.55.3.2 FX3MCSDataDeviceIdOffset const uint32_t FX3MCSDataDeviceIdOffset = 0x4 [static]
- 11.55.3.3 FX3MCSDataIFB1ImageOffset const uint32_t FX3MCSDataIFB1ImageOffset = 0x2C [static]
- 11.55.3.4 FX3MCSDatalFB2ImageOffset const uint32_t FX3MCSDatalFB2ImageOffset = 0xC [static]
- 11.55.3.5 FX3MCSDataVersionOffset const uint32_t FX3MCSDataVersionOffset = 0x8 [static]

11.56 CMcsUsbFunctionNet Class Reference

Inheritance diagram for CMcsUsbFunctionNet:

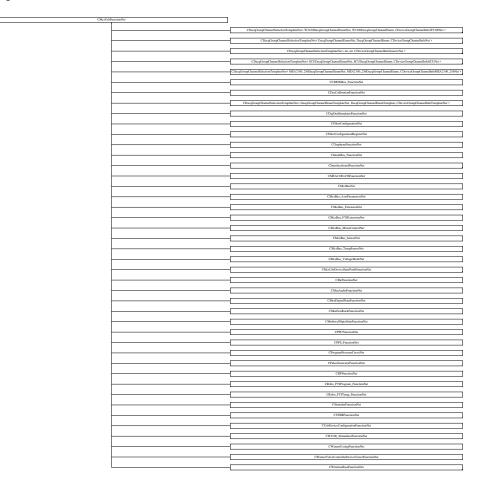

Public Member Functions

- CMcsUsbFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CMcsUsbFunctionNet (void)
- !CMcsUsbFunctionNet ()
- void ThrowCUsbExceptionNetOnError (uint32_t status)

Protected Member Functions

• CMcsUsbFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ mcsusbfunction)

Protected Attributes

- CMcsUsbNet ^ m_pMcsUsb
- CMcsUsbFunction * m_pMcsUsbFunction

11.56.1 Constructor & Destructor Documentation

```
11.56.1.1 CMcsUsbFunctionNet() [1/2] CMcsUsbFunctionNet (
            CMcsUsbNet^ mcsusb )
11.56.1.2 ~CMcsUsbFunctionNet() virtual ~CMcsUsbFunctionNet (
            void ) [virtual]
11.56.1.3 "!CMcsUsbFunctionNet() !CMcsUsbFunctionNet ()
11.56.1.4 CMcsUsbFunctionNet() [2/2] CMcsUsbFunctionNet (
            CMcsUsbNet^ mcsusb,
            CMcsUsbFunctionPointerContainer^{\wedge} mcsusbfunction) [protected]
11.56.2 Member Function Documentation
11.56.2.1 ThrowCUsbExceptionNetOnError() void ThrowCUsbExceptionNetOnError (
            uint32_t status )
11.56.3 Member Data Documentation
11.56.3.1 m_pMcsUsb CMcsUsbNet ^ m_pMcsUsb [protected]
11.56.3.2 m_pMcsUsbFunction CMcsUsbFunction* m_pMcsUsbFunction [protected]
11.57 CMcsUsbFunctionPointerContainer Class Reference
11.58 CMcsUsbListEntryNet Class Reference
```

McsUsbListEntryNet identifies a connected device.

Public Member Functions

- ∼CMcsUsbListEntryNet ()
- virtual bool Equals (Object[^] obj) override

Checks weather two CMcsUsbListEntryNet represent the same USB device.

void SetStringFormat (String ^ format)

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

virtual String \(^{\text{ToString}}\) () override

Static Public Member Functions

static CMcsUsbListEntryNet ^ GetEntry ()

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

• static CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

static CMcsUsbListEntryNet ^ GetEntry (DeviceEnumNet McsUsbDevice, unsigned int index)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

static unsigned int GetEntryCount ()

Returns the number of devices connected to the computer.

static unsigned int GetEntryCount (DeviceEnumNet McsUsbDevice)

Returns the number of devices connected to the computer.

Properties

• String[^] Manufacturer [get]

The Manufacturer ID of the device represented by this CMcsUsbListEntryNet.

String[^] Product [get]

The Product ID of the device represented by this CMcsUsbListEntryNet.

• String DeviceName [get]

The device name of the device represented by this CMcsUsbListEntryNet.

• String[^] SerialNumber [get]

The serial number of the device represented by this CMcsUsbListEntryNet.

• String HwVersion [get]

The hardware revision of the device represented by this CMcsUsbListEntryNet.

DeviceIdNet^ DeviceId [get]

11.58.1 Detailed Description

McsUsbListEntryNet identifies a connected device.

11.58.2 Constructor & Destructor Documentation

11.58.2.1 ~ CMcsUsbListEntryNet() ~ CMcsUsbListEntryNet ()

11.58.3 Member Function Documentation

```
11.58.3.1 Equals() virtual bool Equals (
Object^ obj ) [override], [virtual]
```

Checks weather two CMcsUsbListEntryNet represent the same USB device.

obj The CMcsUsbListEntryNet to compare with.

```
11.58.3.2 GetEntry() [1/3] static CMcsUsbListEntryNet ^ GetEntry ( ) [static]
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

```
11.58.3.3 GetEntry() [2/3] static CMcsUsbListEntryNet ^ GetEntry (

DeviceEnumNet McsUsbDevice ) [static]
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

```
McsUsbDevice | Specifies the type of devices to look for.
```

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

```
11.58.3.4 GetEntry() [3/3] static CMcsUsbListEntryNet ^ GetEntry (

DeviceEnumNet McsUsbDevice,

unsigned int index ) [static]
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

McsUsbDevice | Specifies the type of devices to look for.

index	number of the entry to use.
-------	-----------------------------

Returns

A CMcsUsbListEntryNet to be used to connect to the device.

11.58.3.5 GetEntryCount() [1/2] static unsigned int GetEntryCount () [static]

Returns the number of devices connected to the computer.

Returns

The number of devices.

```
11.58.3.6 GetEntryCount() [2/2] static unsigned int GetEntryCount (

DeviceEnumNet McsUsbDevice ) [static]
```

Returns the number of devices connected to the computer.

Parameters

McsUsbDevice	Specifies the type of devices to look for.
--------------	--

Returns

The number of devices.

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

format A String containing the format template. Default is "%N (%S)".

```
11.58.3.8 ToString() virtual String ^{\wedge} ToString ( ) [override], [virtual]
```

11.58.4 Property Documentation

```
11.58.4.1 DeviceId DeviceIdNet^ DeviceId [get]
```

```
11.58.4.2 DeviceName String^ DeviceName [get]
```

The device name of the device represented by this CMcsUsbListEntryNet.

```
\textbf{11.58.4.3} \quad \textbf{HwVersion} \quad \texttt{String}^{\wedge} \; \texttt{HwVersion} \quad \texttt{[get]}
```

The hardware revision of the device represented by this CMcsUsbListEntryNet.

```
11.58.4.4 Manufacturer String^ Manufacturer [get]
```

The Manufacturer ID of the device represented by this CMcsUsbListEntryNet.

```
11.58.4.5 Product String^ Product [get]
```

The Product ID of the device represented by this ${\sf CMcsUsbListEntryNet}.$

```
11.58.4.6 SerialNumber String SerialNumber [get]
```

The serial number of the device represented by this CMcsUsbListEntryNet.

11.59 CMcsUsbListNet Class Reference

Class to handle a list of connected MCS USB devices.

Public Member Functions

CMcsUsbListNet (DeviceEnumNet McsUsbDevice)

Initializes a new instance of CMcsUsbListNet class.

CMcsUsbListNet (array< DeviceIdNet[^]> DeviceIdList)

Initializes a new instance of CMcsUsbListNet class.

∼CMcsUsbListNet ()

Destructor: called by Dispose()

!CMcsUsbListNet ()

Finalizer: called by GC before collecting

void SetStringFormat (String ^ format)

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

• uint32 t GetNumberOfDevices ()

Gets the number of devices currently in the list.

CMcsUsbListEntryNet ^ GetUsbListEntry (unsigned int index)

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

array< CMcsUsbListEntryNet[^]> [^] GetUsbListEntries ()

Returns all entries from the list of USB Devices connected to the computer.

bool IsDeviceTypeOf (CMcsUsbListEntryNet[^] entry, DeviceEnumNet McsUsbDevice)

Properties

• uint32_t Count [get]

Gets the number of devices currently in the list.

Events

- OnDeviceArrivalRemoval^ DeviceArrival
- OnDeviceArrivalRemoval^ DeviceRemoval

11.59.1 Detailed Description

Class to handle a list of connected MCS USB devices.

11.59.2 Constructor & Destructor Documentation

```
11.59.2.1 CMcsUsbListNet() [1/2] CMcsUsbListNet (

DeviceEnumNet McsUsbDevice )
```

Initializes a new instance of CMcsUsbListNet class.

```
11.59.2.2 CMcsUsbListNet() [2/2] CMcsUsbListNet (

array< DeviceIdNet^>^ DeviceIdList )
```

Initializes a new instance of CMcsUsbListNet class.

```
11.59.2.3 ~CMcsUsbListNet() ~CMcsUsbListNet ()
```

Destructor: called by Dispose()

11.59.2.4 "!CMcsUsbListNet() !CMcsUsbListNet ()

Finalizer: called by GC before collecting

11.59.3 Member Function Documentation

11.59.3.1 GetNumberOfDevices() uint32_t GetNumberOfDevices ()

Gets the number of devices currently in the list.

Returns

The number of devices currently in the list.

```
11.59.3.2 GetUsbListEntries() array<CMcsUsbListEntryNet^> ^ GetUsbListEntries ( )
```

Returns all entries from the list of USB Devices connected to the computer.

```
11.59.3.3 GetUsbListEntry() CMcsUsbListEntryNet ^ GetUsbListEntry ( unsigned int index )
```

Returns one CMcsUsbListEntryNet from the list of USB Devices connected to the computer.

Parameters

index | number of the entry to use.

```
11.59.3.5 SetStringFormat() void SetStringFormat (
String ^ format )
```

Specify the text the CMcsUsbListEntryNet.ToString() function should return. The special code N expands to the device name and S expands to the serial number of the device.

Parameters

	format	A String containing the format template. Default is "%N (%S)".
--	--------	--

11.59.4 Property Documentation

```
11.59.4.1 Count uint32_t Count [get]
```

Gets the number of devices currently in the list.

11.59.5 Event Documentation


11.59.5.1 DeviceArrival OnDeviceArrivalRemoval^ DeviceArrival

11.59.5.2 DeviceRemoval OnDeviceArrivalRemoval^ DeviceRemoval

11.60 CMcsUsbNet Class Reference

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

Inheritance diagram for CMcsUsbNet:

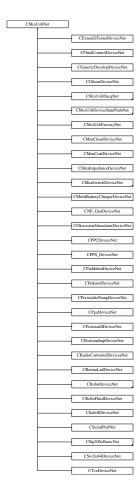

Public Member Functions

• CMcsUsbNet ()

Initializes a new instance of the base class to handle MCS USB devices.

CMcsUsbNet (McsBusTypeEnumNet bustype)

Initializes a new instance of the base class to handle MCS USB devices.

- virtual ∼CMcsUsbNet ()
- !CMcsUsbNet ()
- DeviceEnumNet GetDeviceEnum ()
- virtual uint32_t Connect ()

Opens a connection to the device.

virtual uint32 t Connect (unsigned int LockMask)

Opens a connection to the device.

virtual uint32_t Connect (CMcsUsbListEntryNet[^] entry)

Opens a connection to the device.

• virtual uint32_t Connect (CMcsUsbListEntryNet^ entry, unsigned int LockMask)

Opens a connection to the device.

- virtual uint32_t GetStatus ([System::Runtime::InteropServices::Out]uint32_t% iStatus)
- virtual bool IsConnected ()

Check if a device is Connected.

virtual void Disconnect ()

Disconnect from a device.

- CMcsUsbListEntryNet ^ GetUsbListEntry ()
- virtual String \(^\) GetSerialNumber ()

Query the Serial Number of the device.

- DriverVersionNet [^] GetVersion ()
- DriverVersionNet [^] GetVersion (CFirmwareDestinationNet dest)
- DeviceIdNet ^ GetDeviceId ()
- void MultibootSelectImage (unsigned int sector)

Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.

String \(^\) MultibootGetImageId (unsigned int sector)

Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.

uint32 t MultibootGetCypressImageId (unsigned int sector)

Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.

uint32_t MultibootGetSelectedImage ()

Gets sector index of selected FPGA boot image on IFB

• uint32_t GetMea21UsbPort ()

Gets the USB port if an IFB that is used by this connection

HeadstageIdEnumNet GetHeadstageID (uint32_t headstage)

Gets the ID of a connected headstage.

bool GetHeadstagePresent (uint32 t headstage)

queries whether a headstage is present

bool GetHeadstageActive (uint32 t headstage)

queries whether a headstage is active

• void RescanHeadstage (uint32 t headstage)

rescans and activates a headstage

- array< BYTE > ^ GetSoftwareKey (unsigned int index)
- void SetSoftwareKey (unsigned int index, array< BYTE >^ buffer)
- void RemoveSoftwareKey (unsigned int index)
- void AddSoftwareKey (String[^] key)
- bool EmptyKey (String[^] key)
- bool ValidKey (String[^] key, [System::Runtime::InteropServices::Out]String[^]% serial_number)
- bool HasSoftwareKey (uint8_t ProgrammID, uint8_t majorversion)
- bool HasSoftwareKey (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- String \(^\) GetSoftwareKeyString (uint8 t ProgrammID, uint8 t majorversion)
- String ^ GetSoftwareKeyString (SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID, uint8_t majorversion)
- bool IsDeviceHighSpeedCapable ()
- bool IsDeviceHighSpeed ()
- McsUsbSpeedEnumNet GetDeviceCapableSpeed ()
- McsUsbSpeedEnumNet GetDeviceSpeed ()

Query the Connection Speed of the device.

- unsigned int TxnTestMemoryWrite (unsigned short index)
- unsigned int TxnTestMemoryReadAndCheck (unsigned short index)
- void TxnSetSerialNumber (unsigned int number)
- unsigned int TxnGetSerialNumber ()
- unsigned int ReadRegister (unsigned int reg)
- array< uint32_t > ^ ReadRegister (unsigned int reg, int length)

- unsigned int ReadRegister32 (unsigned int adr)
- unsigned int ReadRegisterTimeSlot (unsigned int reg, int TimeSlot)
- void WriteRegister (unsigned int reg, unsigned int value)
- void WriteRegisterValue (unsigned int reg, unsigned int value)
- void WriteRegister32 (unsigned int adr, unsigned int value)
- void WriteRegister (unsigned int reg, array< unsigned int >^ values)
- void WriteRegisterArray (unsigned int reg, array< unsigned int >^ values)
- · void WriteRegisterTimeSlot (unsigned int reg, unsigned int value, int TimeSlot)
- void WriteRegisterTimeSlot (unsigned int reg, array< unsigned int >^ values, int TimeSlot)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::← InteropServices::Out]uint32_t% DMA_value)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::
 —
 InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize)
- bool ReadEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, [System::Runtime::
 —
 InteropServices::Out]uint32_t% DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void WriteEepromRegisterPreconfig (uint32 t EEPROMBase, uint32 t DMA reg, uint32 t DMA value)
- void WriteEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize)
- void WriteEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t DMA_value, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize)
- void EraseEepromRegisterPreconfig (uint32_t EEPROMBase, uint32_t DMA_reg, uint32_t EEPROMSize, uint32_t EepromStartAddress)
- unsigned int GetLastUSBError ()
- void ThrowCUsbExceptionNetOnError (uint32 t status)
- bool GetDeviceCannotStallOutRequests ()
- String ^ GetHardwareRevision ()
- unsigned int GetFirmwareVersion (CFirmwareDestinationNet destination)

Gets the firmware version for the destination.

- uint8_t GetNumConfigurations ()
- uint8_t GetConfiguration ()
- void SetConfiguration (uint8 t config)
- uint32 t GetDeviceRootHubVendorID ()

Gets the Vendor ID of the USB root hub the device is connected to.

• UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ()

Gets the Vendor ID of the USB root hub the device is connected to.

String ^ GetDeviceRootHubVendorName ()

Gets the Vendor Name of the USB root hub the device is connected to.

void EnableExceptions (bool enable)

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.

- bool IsExceptionsEnabled ()
- uint32_t GetStatusOfLastCommand ()

Gets the status of the last call to the McsUsb Library.

- uint32 t CyclePort ()
- void AssociateToThis (CMcsUsbNet[^] device)

Static Public Member Functions

static String \(^{\text{GetErrorText}}\) (unsigned int Status)

Gets the error text string that belongs to a status number.

Static Public Attributes

```
static const uint32_t Status_Crc = (0xE0100001L)

    static const uint32_t Status_Btstuff = (0xE0100002L)

    static const uint32 t Status DataToggleMismatch = (0xE0100003L)

static const uint32_t Status_Stall = (0xE0100004L)

    static const uint32 t Status DevNotResponding = (0xE0100005L)

    static const uint32 t Status PidCheckFailure = (0xE0100006L)

    static const uint32_t Status_UnexpectedPid = (0xE0100007L)

    static const uint32 t Status DataOverrun = (0xE0100008L)

    static const uint32_t Status_DataUnderrun = (0xE0100009L)

    static const uint32 t Status BufferOverrun = (0xE010000CL)

    static const uint32 t Status BufferUnderrun = (0xE010000DL)

• static const uint32_t Status_NotAccessed = (0xE010000FL)

    static const uint32_t Status_Fifo = (0xE0100010L)

    static const uint32 t Status EndpointHalted = (0xE0100030L)

    static const uint32 t Status NoMemory = (0xE0100100L)

    static const uint32_t Status_InvalidUrbFunction = (0xE0100200L)

    static const uint32_t Status_InvalidParameter = (0xE0100300L)

    static const uint32 t Status InvalidDeviceHandle = (0xE0100013L)

    static const uint32 t Status InvalidHandle = (0xE0100012L)

    static const uint32_t Status_ErrorBusy = (0xE0100400L)

    static const uint32 t Status RequestFailed = (0xE0100500L)

    static const uint32_t Status_InvalidPipeHandle = (0xE0100600L)

    static const uint32_t Status_NoBandwidth = (0xE0100700L)

    static const uint32 t Status InternalHcError = (0xE0100800L)

• static const uint32 t Status ErrorShortTransfer = (0xE0100900L)

    static const uint32 t Status BadStartFrame = (0xE0100A00L)

    static const uint32_t Status_IsochRequestFailed = (0xE0100B00L)

• static const uint32 t Status FrameControlOwned = (0xE0100C00L)

    static const uint32_t Status_ControlNotOwned = (0xE0100D00L)

    static const uint32_t Status_Canceled = (0xE0110000L)

    static const uint32 t Status Canceling = (0xE0120000L)

    static const uint32 t Status AlreadyConfigured = (0xE0110001L)

    static const uint32_t Status_Unconfigured = (0xE0110002L)

    static const uint32_t Status_NoSuchDevice = (0xE01F0002L)

    static const uint32 t Status DeviceNotFound = (0xE01F0003L)

    static const uint32 t Status NotSupported = (0xE01F0005L)

    static const uint32 t Status IoPending = (0xE01F0006L)

    static const uint32_t Status_IoTimeout = (0xE01F0007L)

• static const uint32 t Status DeviceRemoved = (0xE01F0008L)

    static const uint32 t Status PipeNotLinked = (0xE01F0009L)

    static const uint32 t Status ConnectedPipes = (0xE01F000AL)

    static const uint32_t Status_DeviceLocked = (0xE01F0010L)

    static const uint32_t Status_RequestMutexTimeout = (0xE01F0020L)

    static const uint32 t Status RequestMutexFailed = (0xE01F0021L)

    static const uint32_t Status_LastUsbErrorMismatch = (0xE01F0022L)

• static const uint32 t WPAError ScanningIsPending = ( (0xA0220000L) | 0x0036 )
```

Properties

• virtual String Serial Number [get]

11.60.1 Detailed Description

Base class to handle MCS USB devices. All device classes are derived from this class. Functionality that is provided by all MCS devices is handled by this class.

11.60.2 Constructor & Destructor Documentation

```
11.60.2.1 CMcsUsbNet() [1/2] CMcsUsbNet ( )
```

Initializes a new instance of the base class to handle MCS USB devices.

Initializes a new instance of the base class to handle MCS USB devices.

Parameters

bustype	Type of device to use, either USB or PCI.
---------	---

```
11.60.2.3 \sim CMcsUsbNet() virtual \sim CMcsUsbNet () [virtual]
```

```
11.60.2.4 "!CMcsUsbNet() !CMcsUsbNet ()
```

11.60.3 Member Function Documentation

```
11.60.3.1 AddSoftwareKey() void AddSoftwareKey ( String^{\land} key )
```

```
11.60.3.2 AssociateToThis() void AssociateToThis ( CMcsUsbNet^ device )
```

```
11.60.3.3 Connect() [1/4] virtual uint32_t Connect ( ) [virtual]
```

Opens a connection to the device.

Returns

Error Status. 0 on success.

```
11.60.3.4 Connect() [2/4] virtual uint32_t Connect (
CMcsUsbListEntryNet^ entry ) [virtual]
```

Opens a connection to the device.

Parameters

	entry	The Device List Entry for the device to be connected.
--	-------	---

Returns

Error Status. 0 on success.

```
11.60.3.5 Connect() [3/4] virtual uint32_t Connect (

CMcsUsbListEntryNet^ entry,

unsigned int LockMask ) [virtual]
```

Opens a connection to the device.

Parameters

entry	The Device List Entry for the device to be connected.
LockMask	The Lock Mask for this connection.

Returns

Error Status. 0 on success.

Opens a connection to the device.

LockMask	The Lock Mask for this connection.
----------	------------------------------------

Returns

Error Status. 0 on success.

```
 \textbf{11.60.3.7} \quad \textbf{CyclePort()} \quad \texttt{uint32\_t CyclePort ()} \\
```

```
11.60.3.8 Disconnect() virtual void Disconnect ( ) [virtual]
```

Disconnect from a device.

```
11.60.3.9 EmptyKey() bool EmptyKey ( String^{\wedge} key )
```

Enables or Disables Exceptions for calls to McsUsb Devices. If Exceptions are disabled, the return value of a command can be queries with the GetStatusOfLastCommand call instead.

Parameters

enable True to enable Exceptions, False to disable.

```
11.60.3.11 EraseEepromRegisterPreconfig() [1/3] void EraseEepromRegisterPreconfig ( uint32_t EEPROMBase, uint32_t DMA_reg )
```

11.60.3.12 EraseEepromRegisterPreconfig() [2/3] void EraseEepromRegisterPreconfig (

```
uint32_t EEPROMBase,
uint32_t DMA_reg,
uint32_t EEPROMSize )
```

```
11.60.3.13 EraseEepromRegisterPreconfig() [3/3] void EraseEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t EEPROMSize,
             uint32_t EepromStartAddress )
11.60.3.14 GetConfiguration() uint8_t GetConfiguration ( )
11.60.3.15 GetDeviceCannotStallOutRequests() bool GetDeviceCannotStallOutRequests ( )
11.60.3.16 GetDeviceCapableSpeed() McsUsbSpeedEnumNet GetDeviceCapableSpeed ( )
11.60.3.17 GetDeviceEnum() DeviceEnumNet GetDeviceEnum ( )
11.60.3.18 GetDeviceId() DeviceIdNet ^ GetDeviceId ( )
11.60.3.19 GetDeviceRootHubVendorEnum() UsbVendorIdEnumNet GetDeviceRootHubVendorEnum ( )
Gets the Vendor ID of the USB root hub the device is connected to.
Returns
     An enum which enumerates the PCI Vendor ID.
11.60.3.20 GetDeviceRootHubVendorID() uint32_t GetDeviceRootHubVendorID ( )
Gets the Vendor ID of the USB root hub the device is connected to.
Returns
```

The PCI Vendor ID, 0x8086 for Intel, 0x1912 for Renesas, 0x1b21 for ASMedia.

11.60.3.21 GetDeviceRootHubVendorName() String ^ GetDeviceRootHubVendorName ()

Gets the Vendor Name of the USB root hub the device is connected to.

Returns

The PCI Vendor Name, either "Intel", "Renesas", "ASMedia" or "unknown".

11.60.3.22 GetDeviceSpeed() McsUsbSpeedEnumNet GetDeviceSpeed ()

Query the Connection Speed of the device.

Returns

0 for Low-Speed, 1 for Full-Speed, 2 for High-Speed and 3 for SuperSpeed.

11.60.3.23 GetErrorText() static String ^ GetErrorText (unsigned int Status) [static]

Gets the error text string that belongs to a status number.

Parameters

Status	The status number you want the text for.
--------	--

Returns

The error text string that belongs to the status number.

11.60.3.24 **GetFirmwareVersion()** unsigned int GetFirmwareVersion (CFirmwareDestinationNet destination)

Gets the firmware version for the destination.

Parameters

destination	The destination to be queried.

Returns

The firmware version as a 32 bit number, the upper 16 bit contain the majaor version number, the lower 16 bit the minor version number.

11.60.3.25 GetHardwareRevision() String $^{\wedge}$ GetHardwareRevision ()

```
11.60.3.26 GetHeadstageActive() bool GetHeadstageActive ( uint32_t headstage )
```

queries whether a headstage is active

Parameters

in	headstage	the headstage number (0 or 1)
----	-----------	-------------------------------

Returns

true if the headstage is active

11.60.3.27 **GetHeadstageID()** HeadstageIdEnumNet GetHeadstageID (uint32_t headstage)

Gets the ID of a connected headstage.

Parameters

	in	headstage	the headstage number (0 or 1)
--	----	-----------	-------------------------------

Returns

enumerated Headstage ID

11.60.3.28 GetHeadstagePresent() bool GetHeadstagePresent (uint32_t headstage)

queries whether a headstage is present

Parameters

in	headstage	the headstage number (0 or 1)

Returns

true if the headstage is present

```
11.60.3.29 GetLastUSBError() unsigned int GetLastUSBError ()
11.60.3.30 GetMea21UsbPort() uint32_t GetMea21UsbPort ( )
Gets the USB port if an IFB that is used by this connection
Returns
     number of used port; range: 0..1
11.60.3.31 GetNumConfigurations() uint8_t GetNumConfigurations ( )
11.60.3.32 GetSerialNumber() virtual String ^ GetSerialNumber ( ) [virtual]
Query the Serial Number of the device.
Returns
     The Serial Number.
11.60.3.33 GetSoftwareKey() array<BYTE> ^ GetSoftwareKey (
             unsigned int index )
11.60.3.34 GetSoftwareKeyString() [1/2] String ^ GetSoftwareKeyString (
             SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
             uint8_t majorversion )
11.60.3.35 GetSoftwareKeyString() [2/2] String ^{\land} GetSoftwareKeyString (
             uint8_t ProgrammID,
             uint8_t majorversion )
11.60.3.36 GetStatus() virtual uint32_t GetStatus (
             [System::Runtime::InteropServices::Out] uint32_t% iStatus ) [virtual]
```

```
11.60.3.37 GetStatusOfLastCommand() uint32_t GetStatusOfLastCommand ( )
```

Gets the status of the last call to the McsUsb Library.

Returns

The Error Status of the last McsUsb command. 0 on success.

```
11.60.3.38 GetUsbListEntry() CMcsUsbListEntryNet ^ GetUsbListEntry ( )
11.60.3.39 GetVersion() [1/2] DriverVersionNet ^ GetVersion ( )
11.60.3.40 GetVersion() [2/2] DriverVersionNet ^ GetVersion (
             CFirmwareDestinationNet dest )
11.60.3.41 HasSoftwareKey() [1/2] bool HasSoftwareKey (
             SoftwareKeyProgrammIdsNet::ProgrammIdsNet ProgrammID,
             uint8\_t majorversion)
11.60.3.42 HasSoftwareKey() [2/2] bool HasSoftwareKey (
             uint8_t ProgrammID,
             uint8_t majorversion )
11.60.3.43 |sConnected() virtual bool IsConnected ( ) [virtual]
Check if a device is Connected.
Returns
     true if the device is connected.
```

11.60.3.44 IsDeviceHighSpeed() bool IsDeviceHighSpeed ()

```
11.60.3.45 | IsDeviceHighSpeedCapable() bool IsDeviceHighSpeedCapable ( )
11.60.3.46 IsExceptionsEnabled() bool IsExceptionsEnabled ( )
11.60.3.47 MultibootGetCypressImageId() uint32_t MultibootGetCypressImageId (
              unsigned int sector )
Query the multiboot image id of the device located in specified sector (range: 0..9 0..9) of IFB2 Cypress.
Returns
     The magic ident code of the image.
11.60.3.48 MultibootGetImageId() String ^ MultibootGetImageId (
              unsigned int sector )
Query the multiboot image id of the device located in specified sector (range: 0..2 / 0..9) of IFB1 / IFB2 FPGA.
Returns
     The magic ident code of the image.
11.60.3.49 MultibootGetSelectedImage() uint32_t MultibootGetSelectedImage ( )
Gets sector index of selected FPGA boot image on IFB
Returns
     Sector index of image; range: 0..2
11.60.3.50 MultibootSelectImage() void MultibootSelectImage (
              unsigned int sector )
Select the multiboot image specified by "sector" (range: 0..2) for IFB FPGA.
```

Returns

Throws exception on error.

```
\textbf{11.60.3.51} \quad \textbf{ReadEepromRegisterPreconfig() [1/3]} \quad \texttt{bool ReadEepromRegisterPreconfig ()} \quad \textbf{[1/3]} \quad \textbf{
                                                 uint32_t EEPROMBase,
                                                 uint32_t DMA_reg,
                                                 [System::Runtime::InteropServices::Out] uint32_t% DMA_value )
11.60.3.52 ReadEepromRegisterPreconfig() [2/3] bool ReadEepromRegisterPreconfig (
                                                 uint32_t EEPROMBase,
                                                 uint32_t DMA_reg,
                                                 [System::Runtime::InteropServices::Out] uint32_t% DMA_value,
                                                 uint32_t EEPROMSize )
11.60.3.53 ReadEepromRegisterPreconfig() [3/3] bool ReadEepromRegisterPreconfig (
                                                uint32_t EEPROMBase,
                                                uint32_t DMA_reg,
                                                 [System::Runtime::InteropServices::Out] uint32_t% DMA_value,
                                                 uint32_t EEPROMSize,
                                                 uint32_t EepromStartAddress )
11.60.3.54 ReadRegister() [1/2] unsigned int ReadRegister (
                                                 unsigned int reg )
11.60.3.55 ReadRegister() [2/2] array<uint32_t> ^ ReadRegister (
                                                 unsigned int reg,
                                                 int length )
11.60.3.56 ReadRegister32() unsigned int ReadRegister32 (
                                                unsigned int adr )
11.60.3.57 ReadRegisterTimeSlot() unsigned int ReadRegisterTimeSlot (
                                                unsigned int reg,
                                                int TimeSlot )
11.60.3.58 RemoveSoftwareKey() void RemoveSoftwareKey (
                                                 unsigned int index)
11.60.3.59 RescanHeadstage() void RescanHeadstage (
                                                 uint32_t headstage )
rescans and activates a headstage
```

in	headstage	the headstage number ((0 or 1)	
----	-----------	------------------------	----------	--

```
11.60.3.60 SetConfiguration() void SetConfiguration (
              uint8_t config )
\textbf{11.60.3.61} \quad \textbf{SetSoftwareKey()} \quad \texttt{void SetSoftwareKey ()}
              unsigned int index,
              array < BYTE >^{\wedge} buffer )
\textbf{11.60.3.62} \quad \textbf{ThrowCUsbExceptionNetOnError()} \quad \texttt{void ThrowCUsbExceptionNetOnError} \quad \textbf{(}
              uint32_t status )
11.60.3.63 TxnGetSerialNumber() unsigned int TxnGetSerialNumber ()
11.60.3.64 TxnSetSerialNumber() void TxnSetSerialNumber (
              unsigned int number )
11.60.3.65 TxnTestMemoryReadAndCheck() unsigned int TxnTestMemoryReadAndCheck (
               unsigned short index )
11.60.3.66 TxnTestMemoryWrite() unsigned int TxnTestMemoryWrite (
               unsigned short index)
11.60.3.67 ValidKey() [1/2] bool ValidKey (
               String^{\wedge} key,
               [System::Runtime::InteropServices::Out] String^{8} serial_number)
```

```
11.60.3.68 ValidKey() [2/2] bool ValidKey (
             String^{\wedge} key,
             uint8_t ProgrammID,
             uint8_t majorversion,
             [System::Runtime::InteropServices::Out] String^{8} serial_number)
11.60.3.69 WriteEepromRegisterPreconfig() [1/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value )
11.60.3.70 WriteEepromRegisterPreconfig() [2/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value,
             uint32_t EEPROMSize )
11.60.3.71 WriteEepromRegisterPreconfig() [3/3] void WriteEepromRegisterPreconfig (
             uint32_t EEPROMBase,
             uint32_t DMA_reg,
             uint32_t DMA_value,
             uint32_t EEPROMSize,
             uint32_t EepromStartAddress )
11.60.3.72 WriteRegister() [1/2] void WriteRegister (
             unsigned int reg,
             array< unsigned int >^{\wedge} values)
11.60.3.73 WriteRegister() [2/2] void WriteRegister (
             unsigned int reg,
             unsigned int value )
11.60.3.74 WriteRegister32() void WriteRegister32 (
             unsigned int adr,
             unsigned int value )
```

```
11.60.3.75 WriteRegisterArray() void WriteRegisterArray (
             unsigned int reg,
             array< unsigned int >^{\wedge} values )
11.60.3.76 WriteRegisterTimeSlot() [1/2] void WriteRegisterTimeSlot (
             unsigned int reg,
             array< unsigned int >^{\wedge} values,
             int TimeSlot )
\textbf{11.60.3.77} \quad \textbf{WriteRegisterTimeSlot() [2/2]} \quad \texttt{void WriteRegisterTimeSlot ()}
             unsigned int reg,
             unsigned int value,
             int TimeSlot )
11.60.3.78 WriteRegisterValue() void WriteRegisterValue (
             unsigned int reg,
             unsigned int value )
11.60.4 Member Data Documentation
11.60.4.1 Status_AlreadyConfigured const uint32_t Status_AlreadyConfigured = (0xE0110001L)
[static]
11.60.4.2 Status_BadStartFrame const uint32_t Status_BadStartFrame = (0xE0100A00L) [static]
11.60.4.3 Status_Btstuff const uint32_t Status_Btstuff = (0xE0100002L) [static]
11.60.4.4 Status_BufferOverrun const uint32_t Status_BufferOverrun = (0xE010000CL) [static]
11.60.4.5 Status_BufferUnderrun const uint32_t Status_BufferUnderrun = (0xE010000DL) [static]
```

```
11.60.4.6 Status_Canceled const uint32_t Status_Canceled = (0xE0110000L) [static]
11.60.4.7 Status_Canceling const uint32_t Status_Canceling = (0xE0120000L) [static]
11.60.4.8 Status_ConnectedPipes const uint32_t Status_ConnectedPipes = (0xE01F000AL) [static]
11.60.4.9 Status_ControlNotOwned const uint32_t Status_ControlNotOwned = (0xE0100D00L) [static]
11.60.4.10 Status_Crc const uint32_t Status_Crc = (0xE0100001L) [static]
11.60.4.11 Status_DataOverrun const uint32_t Status_DataOverrun = (0xE0100008L) [static]
11.60.4.12 Status_DataToggleMismatch const uint32_t Status_DataToggleMismatch = (0xE0100003L)
[static]
11.60.4.13 Status DataUnderrun const uint32_t Status_DataUnderrun = (0xE0100009L) [static]
11.60.4.14 Status_DeviceLocked const uint32_t Status_DeviceLocked = (0xE01F0010L) [static]
11.60.4.15 Status_DeviceNotFound const uint32_t Status_DeviceNotFound = (0xE01F0003L) [static]
11.60.4.16 Status_DeviceRemoved const uint32_t Status_DeviceRemoved = (0xE01F0008L) [static]
```

```
11.60.4.17 Status_DevNotResponding const uint32_t Status_DevNotResponding = (0xE0100005L)
[static]
11.60.4.18 Status_EndpointHalted const uint32_t Status_EndpointHalted = (0xE0100030L) [static]
11.60.4.19 Status_ErrorBusy const uint32_t Status_ErrorBusy = (0xE0100400L) [static]
11.60.4.20 Status_ErrorShortTransfer const uint32_t Status_ErrorShortTransfer = (0xE0100900L)
[static]
11.60.4.21 Status_Fifo const uint32_t Status_Fifo = (0xE0100010L) [static]
11.60.4.22 Status_FrameControlOwned const uint32_t Status_FrameControlOwned = (0xE0100C00L)
[static]
11.60.4.23 Status_InternalHcError const uint32_t Status_InternalHcError = (0xE0100800L) [static]
11.60.4.24 Status_InvalidDeviceHandle const uint32_t Status_InvalidDeviceHandle = (0xE0100013L)
[static]
11.60.4.25 Status_InvalidHandle const uint32_t Status_InvalidHandle = (0xE0100012L) [static]
11.60.4.26 Status_InvalidParameter const uint32_t Status_InvalidParameter = (0xE0100300L) [static]
```


```
11.60.4.27 Status_InvalidPipeHandle const uint32_t Status_InvalidPipeHandle = (0xE0100600L)
[static]
11.60.4.28 Status_InvalidUrbFunction const uint32_t Status_InvalidUrbFunction = (0xE0100200L)
[static]
11.60.4.29 Status_loPending const uint32_t Status_IoPending = (0xE01F0006L) [static]
11.60.4.30 Status_loTimeout const uint32_t Status_IoTimeout = (0xE01F0007L) [static]
11.60.4.31 Status_IsochRequestFailed const uint32_t Status_IsochRequestFailed = (0xE0100B00L)
[static]
11.60.4.32 Status LastUsbErrorMismatch const uint32_t Status_LastUsbErrorMismatch = (0xE01↔
F0022L) [static]
11.60.4.33 Status_NoBandwidth const uint32_t Status_NoBandwidth = (0xE0100700L) [static]
11.60.4.34 Status_NoMemory const uint32_t Status_NoMemory = (0xE0100100L) [static]
11.60.4.35 Status_NoSuchDevice const uint32_t Status_NoSuchDevice = (0xE01F0002L) [static]
11.60.4.36 Status_NotAccessed const uint32_t Status_NotAccessed = (0xE010000FL) [static]
11.60.4.37 Status_NotSupported const uint32_t Status_NotSupported = (0xE01F0005L) [static]
```

```
11.60.4.38 Status_PidCheckFailure const uint32_t Status_PidCheckFailure = (0xE0100006L) [static]
11.60.4.39 Status_PipeNotLinked const uint32_t Status_PipeNotLinked = (0xE01F0009L) [static]
11.60.4.40 Status_RequestFailed const uint32_t Status_RequestFailed = (0xE0100500L) [static]
11.60.4.41 Status_RequestMutexFailed const uint32_t Status_RequestMutexFailed = (0xE01F0021L)
 [static]
\textbf{11.60.4.42} \quad \textbf{Status\_RequestMutexTimeout} \quad \texttt{const uint} \\ \textbf{32\_t Status\_RequestMutexTimeout} = (0 \\ \texttt{xE01} \\ \leftarrow \\ \textbf{0xE01} \\ \leftarrow \\ \textbf{0x
F0020L) [static]
11.60.4.43 Status_Stall const uint32_t Status_Stall = (0xE0100004L) [static]
11.60.4.44 Status_Unconfigured const uint32_t Status_Unconfigured = (0xE0110002L) [static]
11.60.4.45 Status_UnexpectedPid const uint32_t Status_UnexpectedPid = (0xE0100007L) [static]
11.60.4.46 WPAError_ScanningIsPending const uint32_t WPAError_ScanningIsPending = ( (0x↔
A0220000L) | 0x0036 ) [static]
11.60.5 Property Documentation
11.60.5.1 SerialNumber virtual String^ SerialNumber [get]
```

11.61 CMcsUsbPointerContainer Class Reference

11.62 CMEA2100_256DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CMEA2100_256DacqGroupChannelSelectionNet:

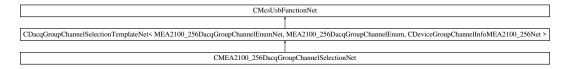

Public Member Functions

CMEA2100_256DacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members

11.62.1 Constructor & Destructor Documentation


```
11.62.1.1 CMEA2100_256DacqGroupChannelSelectionNet() CMEA2100_256DacqGroupChannelSelectionNet (

CMcsUsbNet^ mcsusb )
```

11.63 CMEA2100x256FunctionNet Class Reference

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

Inheritance diagram for CMEA2100x256FunctionNet:

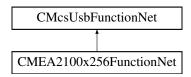

Public Member Functions

CMEA2100x256FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pMEA2100x256←
 FunctionPointerContainer)

Initializes a new instance of the CMEA2100x256FunctionNet class.

- CMEA2100x256FunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CMEA2100x256FunctionNet ()
- !CMEA2100x256FunctionNet ()
- StimulationLayoutConfigurationEnumNet GetLayoutConfiguration ()

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulation

void SetLayoutConfiguration (StimulationLayoutConfigurationEnumNet LayoutConfiguration)

Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulation and the stimulation is the stimulation of the stimulation of the stimulation is the stimulation of the stimulation of the stimulation is the stimulation of the stimul

Additional Inherited Members

11.63.1 Detailed Description

CMEA2100x256FunctionNet is the class to control the MEA2100-256 device needs #include "Stg200xNet.h" to resolve documentation reference

11.63.2 Constructor & Destructor Documentation

```
11.63.2.1 CMEA2100x256FunctionNet() [1/2] CMEA2100x256FunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pMEA2100x256FunctionPointerContainer)
```

Initializes a new instance of the CMEA2100x256FunctionNet class.

```
11.63.2.2 CMEA2100x256FunctionNet() [2/2] CMEA2100x256FunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.63.2.3 ~CMEA2100x256FunctionNet() virtual ~CMEA2100x256FunctionNet ( ) [virtual]
```

```
11.63.2.4 "!CMEA2100x256FunctionNet() !CMEA2100x256FunctionNet ()
```

11.63.3 Member Function Documentation

```
11.63.3.1 GetLayoutConfiguration() StimulationLayoutConfigurationEnumNet GetLayoutConfiguration ()
```

Gets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

Returns

The currently active stimulation layout configuration.


```
11.63.3.2 SetLayoutConfiguration() void SetLayoutConfiguration (
StimulationLayoutConfigurationEnumNet LayoutConfiguration)
```

Sets the stimulation layout configuration. Can be single well, 6-well or 9-well. The number of DAC channels available per well is Mcs::Usb::CStg200xBasicNet::GetNumberOfAnalogChannels divided by Mcs::Usb::CStg200xBasicNet::GetNumberOfStimulationSourcesPerElectrode.

LayoutConfiguration	The new stimulation layout configuration.	
---------------------	---	--

11.64 CMeaAudioFunctionNet Class Reference

Inheritance diagram for CMeaAudioFunctionNet:

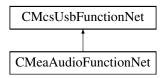

Classes

· struct s_setaudionet

Public Member Functions

- CMeaAudioFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] meaAudioFunction←
 PointerContainer)
- CMeaAudioFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual uint32_t GetNumberOfAudioChannels ()

Gets the number of available audio channels.

virtual uint32_t SetAudioChannels (array< s_setaudionet^>^ channels)

Sets the electrode to monitor and amplification for the audio channels.

- virtual uint32_t SetAudioChannels (array < s_setaudionet^>^ channels, unsigned int virtualDevice)
 Sets the electrode to monitor and amplification for the audio channels.
- virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out]array< s_setaudionet^>^% channels)

Gets the electrode to monitor and amplification for the audio channels.

virtual uint32_t GetAudioChannels ([System::Runtime::InteropServices::Out]array< s_setaudionet^>^% channels, unsigned int virtualDevice)

Gets the electrode to monitor and amplification for the audio channels.

Additional Inherited Members

11.64.1 Constructor & Destructor Documentation

```
11.64.1.1 CMeaAudioFunctionNet() [1/2] CMeaAudioFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ meaAudioFunctionPointerContainer)
```

```
11.64.1.2 CMeaAudioFunctionNet() [2/2] CMeaAudioFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.64.2 Member Function Documentation

Gets the electrode to monitor and amplification for the audio channels.

Parameters

channels Struct which contains the electrode (channel) and amplification on return.

Returns

Error Status. 0 on success.

Gets the electrode to monitor and amplification for the audio channels.

Parameters

channels Struct which contains the electrode (channel) and amplification on return.

Parameters

Returns

Error Status. 0 on success.

11.64.2.3 GetNumberOfAudioChannels() virtual uint32_t GetNumberOfAudioChannels () [virtual]

Gets the number of available audio channels.

Returns

The number of audio channels available, 0 when there are none.

```
11.64.2.4 SetAudioChannels() [1/2] virtual uint32_t SetAudioChannels (

array< s_setaudionet^>^ channels ) [virtual]
```

Sets the electrode to monitor and amplification for the audio channels.

Parameters

Returns

Error Status. 0 on success.

Sets the electrode to monitor and amplification for the audio channels.

Parameters

channels Struct which defines the electrode (channel) and amplific
--

Parameters


Returns

Error Status. 0 on success.

11.65 CMeaCleanDeviceNet Class Reference

CMeaCleanDeviceNet is the class to access the MEA Clean device.

Inheritance diagram for CMeaCleanDeviceNet:

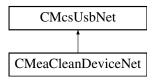

Public Member Functions

• CMeaCleanDeviceNet ()

Initializes a new instance of the CMeaCleanDeviceNet class.

- virtual ∼CMeaCleanDeviceNet ()
- !CMeaCleanDeviceNet ()
- void Start ()

Starts a MEA Clean run.

• void Stop ()

Stops a MEA Clean run.

void SetSlope (uint32_t voltageSlope)

Sets the voltage slope.

void SetCycles (uint32_t cycles)

Sets the number of cycles.

• void SetMinVoltage (int32_t voltageMin)

Sets the lower voltage level.

void SetMaxVoltage (int32_t voltageMax)

Sets the upper voltage level.

• bool IsRunning ()

Gets if the MEA Clean device is running.

• uint32 t GetSlope ()

Gets the voltage slope.

uint32_t GetCycles ()

Gets the number of cycles.

int32_t GetMinVoltage ()

Gets the lower voltage level.

int32_t GetMaxVoltage ()

Gets the upper voltage level

• int32_t GetOutputVoltage ()

Gets the output voltage.

• int32_t GetCycle ()

Gets the current cycle.

Additional Inherited Members

11.65.1 Detailed Description

CMeaCleanDeviceNet is the class to access the MEA Clean device.

11.65.2 Constructor & Destructor Documentation

Initializes a new instance of the CMeaCleanDeviceNet class.

```
11.65.2.2 ~CMeaCleanDeviceNet() virtual ~CMeaCleanDeviceNet () [virtual]
```

```
11.65.2.3 "!CMeaCleanDeviceNet() !CMeaCleanDeviceNet ()
```

11.65.3 Member Function Documentation

```
11.65.3.1 GetCycle() int32_t GetCycle ()
```

Gets the current cycle.

Returns

The cycle number.

```
11.65.3.2 GetCycles() uint32_t GetCycles ()
```

Gets the number of cycles.

Returns

The number of cycles to run for.

```
11.65.3.3 GetMaxVoltage() int32_t GetMaxVoltage ( )
Gets the upper voltage level
Returns
     The upper voltage level in mV.
11.65.3.4 GetMinVoltage() int32_t GetMinVoltage ( )
Gets the lower voltage level.
Returns
     The lower voltage level in mV.
11.65.3.5 GetOutputVoltage() int32_t GetOutputVoltage ()
Gets the output voltage.
Returns
     The output voltage in mV.
11.65.3.6 GetSlope() uint32_t GetSlope ()
Gets the voltage slope.
Returns
     The voltage slope in mV/s.
11.65.3.7 IsRunning() bool IsRunning ( )
Gets if the MEA Clean device is running.
Returns
     "true" when a run is in progress, otherwise "false".
```

Sets the number of cycles.

Sets the upper voltage level.

Parameters

	voltageMax	The upper voltage level in mV (-1.6 1.6 V).
--	------------	---

Sets the lower voltage level.

Parameters

	voltageMin	The lower voltage level in mV (-1.6 1.6 V).
--	------------	---

Sets the voltage slope.

Parameters

voltageSlope The voltage slope in mV/s (range 0	60 V/s).
---	----------

```
11.65.3.12 Start() void Start ( )
```

Starts a MEA Clean run.


11.65.3.13 Stop() void Stop ()

Stops a MEA Clean run.

11.66 CMeaCoatDeviceNet Class Reference

CMeaCoatDeviceNet is the class to access the MEA Coat device.

Inheritance diagram for CMeaCoatDeviceNet:

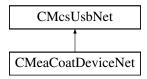

Public Member Functions

CMeaCoatDeviceNet ()

Initializes a new instance of the CMeaCoatDeviceNet class.

- virtual ∼CMeaCoatDeviceNet ()
- !CMeaCoatDeviceNet ()
- void Start ()

Starts a MEA Coat run.

• void Stop ()

Stops a MEA Coat run.

void SetSlope (int32_t currentSlope)

Sets the current slope.

void SetDuration (uint32_t duration)

Sets the duration of a MEA Coat run.

void SetMaxCurrent (uint32_t currentMax)

Sets the limit of the current ramp (absolute value).

void SetOffsetCurrent (int32_t currentOffset)

Sets the offset of the current.

• bool IsRunning ()

Gets if the MEA Clean device is running.

• int32_t GetSlope ()

Gets the current slope.

• uint32_t GetDuration ()

Gets the duration of a MEA Coat run.

• uint32_t GetMaxCurrent ()

Gets the limit of the current ramp (absolute value).

• int32 t GetOffsetCurrent ()

Gets the offset of the current.

• int32_t GetOutputCurrent ()

Gets the output current.

• int32_t GetTimeInPlateau ()

Gets the time in the plateau.

void SetPauseDuration (uint32 t pauseDuration)

Sets the duration of the pause between MEA Coat pulses.

· uint32_t GetPauseDuration ()

Gets the duration of the pause between MEA Coat pulses.

int32_t GetTimeInPause ()

Gets the time in the pause.

```
    void SetCycles (uint32_t cycles)

         Sets the number of cycles.
    • uint32_t GetCycles ()
         Gets the number of cycles.

    int32_t GetCurrentCycle ()

         Gets the current cycle.
Additional Inherited Members
11.66.1 Detailed Description
CMeaCoatDeviceNet is the class to access the MEA Coat device.
11.66.2 Constructor & Destructor Documentation
11.66.2.1 CMeaCoatDeviceNet() CMeaCoatDeviceNet ()
Initializes a new instance of the CMeaCoatDeviceNet class.
11.66.2.2 ~ CMeaCoatDeviceNet() virtual ~ CMeaCoatDeviceNet () [virtual]
11.66.2.3 "!CMeaCoatDeviceNet() !CMeaCoatDeviceNet ()
11.66.3 Member Function Documentation
11.66.3.1 GetCurrentCycle() int32_t GetCurrentCycle ( )
Gets the current cycle.
Returns
     The cycle number.
```

The output current in pA.

```
11.66.3.2 GetCycles() uint32_t GetCycles ()
Gets the number of cycles.
Returns
     The number of cycles to run for.
11.66.3.3 GetDuration() uint32_t GetDuration ( )
Gets the duration of a MEA Coat run.
Returns
     The duration in ms.
11.66.3.4 GetMaxCurrent() uint32_t GetMaxCurrent ( )
Gets the limit of the current ramp (absolute value).
Returns
     The limit of the current ramp in pA (absolute value).
11.66.3.5 GetOffsetCurrent() int32_t GetOffsetCurrent ( )
Gets the offset of the current.
Returns
     The offset of the current in pA.
11.66.3.6 GetOutputCurrent() int32_t GetOutputCurrent ( )
Gets the output current.
Returns
```

11.66.3.7 GetPauseDuration() uint32_t GetPauseDuration ()

Gets the duration of the pause between MEA Coat pulses.

Returns

The duration in ms.

```
11.66.3.8 GetSlope() int32_t GetSlope ()
```

Gets the current slope.

Returns

The current slope in pA/s.

11.66.3.9 GetTimeInPause() int32_t GetTimeInPause ()

Gets the time in the pause.

Returns

The time in the pause in ms.

11.66.3.10 GetTimeInPlateau() int32_t GetTimeInPlateau ()

Gets the time in the plateau.

Returns

The time in the plateau in ms.

```
11.66.3.11 IsRunning() bool IsRunning ()
```

Gets if the MEA Clean device is running.

Returns

"true" when a run is in progress, otherwise "false".

Sets the number of cycles.

cycles The number of cycles to run for (0	9).
---	-----

11.66.3.13 SetDuration() void SetDuration (uint32_t duration)

Sets the duration of a MEA Coat run.

Parameters

tion in ms (range 0 65 s)	The duration in ms	duration
---------------------------	--------------------	----------

11.66.3.14 SetMaxCurrent() void SetMaxCurrent (uint32_t currentMax)

Sets the limit of the current ramp (absolute value).

Parameters

currentMax	The limit of the current ramp in pA (absolute value, 0 18 nA).

Sets the offset of the current.

Parameters

currentOffset	The offset of the current in pA (-10 10 nA).
---------------	--

11.66.3.16 SetPauseDuration() void SetPauseDuration (uint32_t pauseDuration)

Sets the duration of the pause between MEA Coat pulses.

Parameters

nauseDuration	The duration in ms (range	0 65 s)

Sets the current slope.

Parameters

currentSlope	The current slope in pA/s (range -65 65 nA/s).
--------------	--

```
11.66.3.18 Start() void Start ()
```

Starts a MEA Coat run.


```
11.66.3.19 Stop() void Stop ()
```

Stops a MEA Coat run.

11.67 CMeaDeviceNet Class Reference

Base class for MEA data acquisition devices.

Inheritance diagram for CMeaDeviceNet:

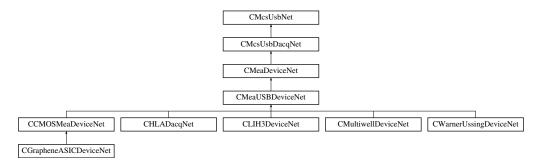

Public Member Functions

- CMeaDeviceNet (McsBusTypeEnumNet bustype)
 - Initializes a new instance of CMeaDeviceNet class.
- CMeaDeviceNet (McsBusTypeEnumNet bustype, OnChannelData[^] channelData, OnError[^] error)
 Initializes a new instance of CMeaDeviceNet class.
- ∼CMeaDeviceNet ()
- virtual int32 t GetGain ()

Gets the amplifier gain of the device.

- int32_t GetEnumerationSpeed ()
- virtual int32_t GetAnalogGain ()

Gets the gain of the analog inputs of the device.

• virtual uint32_t EnableDigitalIn (bool enable, unsigned int virtualDevice)

Enable the digital data word in the datastream.

virtual uint32_t EnableDigitalIn (DigitalDatastreamEnableEnumNet enable, unsigned int virtualDevice)

Enable digital data words in the datastream.

virtual uint32 t EnableTimestamp (bool enable, unsigned int virtualDevice)

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

virtual uint32_t EnableChecksum (bool enable, unsigned int virtualDevice)

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

virtual void SetDigitalOut (unsigned int digout value, int pulselength)

Generate a pulse on the digital output.

• virtual uint32 t SetNumberOfChannels (int NumberOfChannels)

Sets the number of analog channels in the datastream.

virtual uint32 t SetNumberOfChannels (int NumberOfChannels, unsigned int virtualDevice)

Sets the number of analog channels in the datastream.

virtual uint32_t SetNumberOfAnalogChannels (unsigned int NumberOfChannels_HS1, unsigned int NumberOfChannels_HS2, unsigned int NumberOfChannels_DSP, unsigned int NumberOfChannels_IF, unsigned int virtualDevice)

Sets the number of analog channels in the datastream for the MEA2100 device.

virtual uint32_t SetTriggerPeriod (int samples, unsigned int virtualDevice)

Sets the maximum number of samples per trigger.

virtual uint32_t SetTriggerMaskValue (unsigned int mask, unsigned int value, unsigned int virtualDevice)

Defines a pattern on the digital dataword which will start a trigger when found.

Properties

- CMeFunctionNet[^] MeFunctionNet [get]
- CWClassicFunctionNet[^] WClassicFunctionNet [get]
- CW2100 FunctionNet[^] W2100 FunctionNet [get]
- CMeaAudioFunctionNet[^] MeaAudioFunctionNet [get]
- CMeaDigitalDataFunctionNet[^] MeaDigitalDataFunctionNet [get]
- CMeaFeedbackFunctionNet[^] MeaFeedbackFunctionNet [get]
- virtual int Gain [get]

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

virtual int AnalogGain [get]

The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Additional Inherited Members

11.67.1 Detailed Description

Base class for MEA data acquisition devices.

There are two different device types for MEA data aquistion devices. There are the USB-MEA devices and the MC← _Card. In .NET both classes can be accessed by the contructor of the base class CMeaDeviceNet, which contructs the correct underlying C++ class for the USB-MEA device on the one hand or the MC_Card device on the other hand. Through this interface both device types USB-MEA devices and MC_Card devices can be accessed

11.67.2 Constructor & Destructor Documentation

11.67.2.1 CMeaDeviceNet() [1/2] CMeaDeviceNet (McsBusTypeEnumNet bustype)

Initializes a new instance of CMeaDeviceNet class.

Parameters

bustype	Type of device to use, either USB or PCI.
---------	---

11.67.2.2 CMeaDeviceNet() [2/2] CMeaDeviceNet (

```
\begin{tabular}{ll} McsBusTypeEnumNet & bustype, \\ OnChannelData^ & channelData, \\ OnError^ & error \end{tabular} \label{table}
```

Initializes a new instance of CMeaDeviceNet class.

Parameters

bustype	Type of device to use, either USB or PCI.
---------	---

Parameters

channelData	Callback to call when new data is available.
-------------	--

Parameters

error Callback to call when an error occurred.

11.67.2.3 ~CMeaDeviceNet() ~CMeaDeviceNet ()

11.67.3 Member Function Documentation

Enable the checksum data word in the datastream. The checksum is a 32 bit counter and 2x16 bit magic numbers.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable the digital data word in the datastream.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable digital data words in the datastream.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Enable the timestamp data word in the datastream. The timestamp is a 64 bit counter.

Parameters

enable	True to enable, False to disable.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

11.67.3.5 GetAnalogGain() virtual int32_t GetAnalogGain () [virtual]

Gets the gain of the analog inputs of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

```
11.67.3.6 GetEnumerationSpeed() int32_t GetEnumerationSpeed ( )
```

```
11.67.3.7 GetGain() virtual int32_t GetGain ( ) [virtual]
```

Gets the amplifier gain of the device.

Returns

Gain times 1000, a value of 1000 corresponds to a gain of 1.0.

Generate a pulse on the digital output.

out_value Bitmask to set on the digital out.
--

Parameters

pulselength Pulselength in ms

Sets the number of analog channels in the datastream for the MEA2100 device.

Parameters

NumberOfChannels HS1	Number of analog channels from the Headstage 1.	
	,	

Parameters

NumberOfChannels_HS2	Number of analog channels from the Headstage 2.
----------------------	---

Parameters

NumberOfChannels DSP Number of data words from the DSP.

Parameters

NumberOfChannels↔	Number of analog channels from the Interfaceboard.
_IF	

Returns

Error Status. 0 on success.

```
11.67.3.10 SetNumberOfChannels() [1/2] virtual uint32_t SetNumberOfChannels ( int NumberOfChannels ) [virtual]
```

Sets the number of analog channels in the datastream.

Parameters

NumberOfChannels	Number of analog channels.
------------------	----------------------------

Returns

Error Status. 0 on success.

Sets the number of analog channels in the datastream.

Parameters

NumberOfChannels	Number of analog channels.
virtualDevice	virtual device to use.

Returns

Error Status. 0 on success.

Defines a pattern on the digital dataword which will start a trigger when found.

Parameters

	mask	Bits in the digital dataword which are monitored for a match with value.	
--	------	--	--

Parameters

Returns

Error Status. 0 on success.

Sets the maximum number of samples per trigger.

Parameters

Returns

Error Status. 0 on success.

11.67.4 Property Documentation

11.67.4.1 AnalogGain virtual int AnalogGain [get]

The gain of the analog inputs of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

```
11.67.4.2 Gain virtual int Gain [get]
```

The amplifier gain of the device. Value is gain times 1000, a value of 1000 corresponds to a gain of 1.0.

11.67.4.3 MeaAudioFunctionNet CMeaAudioFunctionNet^ MeaAudioFunctionNet [get]

11.67.4.4 MeaDigitalDataFunctionNet CMeaDigitalDataFunctionNet^ MeaDigitalDataFunctionNet [get]

11.67.4.5 MeaFeedbackFunctionNet CMeaFeedbackFunctionNet^ MeaFeedbackFunctionNet [get]


```
11.67.4.6 MeFunctionNet CMeFunctionNet^ MeFunctionNet [get]
```

```
11.67.4.7 W2100_FunctionNet CW2100_FunctionNet^ W2100_FunctionNet [get]
```

11.67.4.8 WClassicFunctionNet CWClassicFunctionNet^ WClassicFunctionNet [get]

11.68 CMeaDigitalDataFunctionNet Class Reference

Inheritance diagram for CMeaDigitalDataFunctionNet:

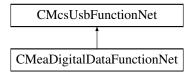

Public Member Functions

- CMeaDigitalDataFunctionNet (CMcsUsbNet[^] mcsusb)
- void SetDigitalData (unsigned int digital_value, unsigned int digital_value_mask)

Generate a value on the digital output.

void SetDigitalData (unsigned int bit_number, bool value)

Generate a value on the digital output.

unsigned int GetDigitalData ()

Get the value of the digital output.

Additional Inherited Members

11.68.1 Constructor & Destructor Documentation

11.68.2 Member Function Documentation

```
11.68.2.1 GetDigitalData() unsigned int GetDigitalData ( )
```

Get the value of the digital output.

Returns

Value on the digital data register.

Generate a value on the digital output.

Parameters

bit_number	Bit number to change.
------------	-----------------------

Parameters

value Bit value.

```
11.68.2.3 SetDigitalData() [2/2] void SetDigitalData (
unsigned int digital_value,
unsigned int digital_value_mask)
```

Generate a value on the digital output.

Parameters


digital_value	Value to set.
---------------	---------------

Parameters

digital_value_ı	<i>mask</i> Mas	sk for change.
-----------------	-------------------	----------------

11.69 CMeaFeedbackFunctionNet Class Reference

Inheritance diagram for CMeaFeedbackFunctionNet:

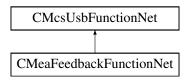

Public Member Functions

- CMeaFeedbackFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] meaFeedback
 FunctionNet)
- CMeaFeedbackFunctionNet (CMcsUsbNet[^] mcsusb)
- void FeedbackSetFeedback (unsigned char on, unsigned short digoutmask, unsigned short diginmask)
- unsigned int FeedbackGetSampleTimerCount ([System::Runtime::InteropServices::Out]unsigned int% CurrentCount, [System::Runtime::InteropServices::Out]unsigned int% LastKnownCount, [System::Runtime ::InteropServices::Out]bool% On)
- void FeedbackSetDigitalMapping (unsigned short channel, unsigned short outmapping, unsigned short inmapping)
- void FeedbackSetFilterParameter (unsigned char filter, array< short >^ parameters)
- void FeedbackSetFilterParameter32 (unsigned char filter, array< int >^ parameters)
- void FeedbackSetIIRFilterParameter (unsigned char filter, int length, array< double >^ parameters)
- void FeedbackSetMkFilter (unsigned char filter, String^ filtertype, double cheb_ribble, String^ passtype, int order, double alpha1, double alpha2)
- void FeedbackSetCreateFilter (unsigned char filter, CCreateFilterNet[^] createFilter)
- void FeedbackSetChannelFilter (short channel, char filter)

- void FeedbackSetGlobalChannelFilter (char filter, unsigned short firstchannel, unsigned short lastchannel)
- void FeedbackSetFilterOff ()
- void FeedbackSetNumberOfSpikeDetectors (unsigned short number)
- void FeedbackSetSpikeDetectorThreshold (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, int threshold1, int threshold2, short slope)
- void FeedbackSetNumberOfRateCounter (unsigned short number)
- void FeedbackSetRateCounter (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel)
- void FeedbackSetNumberOfRateDetectors (unsigned short number)
- void FeedbackSetRateDetector (unsigned short position, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit, unsigned short pulses, unsigned int duration1, unsigned int duration2)
- void FeedbackSetNumberOfLogics (unsigned short number)
- void FeedbackSetLogic (unsigned short position, array< unsigned short >^ sourcechannel, unsigned short resultchannel, unsigned int lookup)
- void FeedbackSetNumberOfTriggers (unsigned short number)
- void FeedbackSetTrigger (unsigned short position, unsigned short sourcechannel, unsigned short resultchannel, unsigned short trigger, unsigned short totzeit)
- void FeedbackSetAnalogSource (AnalogSourceEnumNet AnalogSource, unsigned int Channels, unsigned int Offset)
- void FeedbackSetTriggerTotzeitFactor (unsigned short trigger_factor, unsigned short totzeit_factor)
- void FeedbackGetTriggerTotzeitFactor ([System::Runtime::InteropServices::Out]unsigned short% trigger_
 factor, [System::Runtime::InteropServices::Out]unsigned short% totzeit factor)

Additional Inherited Members

11.69.1 Constructor & Destructor Documentation

```
11.69.1.1 CMeaFeedbackFunctionNet() [1/2] CMeaFeedbackFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ meaFeedbackFunctionNet )
```

```
11.69.1.2 CMeaFeedbackFunctionNet() [2/2] CMeaFeedbackFunctionNet ( CMcsUsbNet^ mcsusb )
```

11.69.2 Member Function Documentation

```
11.69.2.2 FeedbackGetTriggerTotzeitFactor() void FeedbackGetTriggerTotzeitFactor (
               [System::Runtime::InteropServices::Out] unsigned short% trigger_factor,
               [{\tt System::Runtime::InteropServices::Out}] \ \ unsigned \ \ short \% \ \ totzeit\_factor \ )
11.69.2.3 FeedbackSetAnalogSource() void FeedbackSetAnalogSource (
              AnalogSourceEnumNet AnalogSource,
              unsigned int Channels,
              unsigned int Offset )
11.69.2.4 FeedbackSetChannelFilter() void FeedbackSetChannelFilter (
              short channel,
              char filter )
11.69.2.5 FeedbackSetCreateFilter() void FeedbackSetCreateFilter (
              unsigned char filter,
              CCreateFilterNet^ createFilter )
\textbf{11.69.2.6} \quad \textbf{FeedbackSetDigitalMapping()} \quad \texttt{void FeedbackSetDigitalMapping ()}
              unsigned short channel,
              unsigned short outmapping,
              unsigned short inmapping )
\textbf{11.69.2.7} \quad \textbf{FeedbackSetFeedback()} \quad \texttt{void FeedbackSetFeedback} \quad \textbf{(}
              unsigned char on,
              unsigned short digoutmask,
              unsigned short diginmask )
11.69.2.8 FeedbackSetFilterOff() void FeedbackSetFilterOff ()
11.69.2.9 FeedbackSetFilterParameter() void FeedbackSetFilterParameter (
              unsigned char filter,
              array < short >^{\land} parameters)
```


```
11.69.2.10 FeedbackSetFilterParameter32() void FeedbackSetFilterParameter32 (
               unsigned char filter,
               array< int >^{\wedge} parameters )
\textbf{11.69.2.11} \quad \textbf{FeedbackSetGlobalChannelFilter()} \quad \texttt{void FeedbackSetGlobalChannelFilter} \quad \textbf{(}
               char filter,
               unsigned short firstchannel,
               unsigned short lastchannel )
\textbf{11.69.2.12} \quad \textbf{FeedbackSetIIRFilterParameter()} \quad \texttt{void FeedbackSetIIRFilterParameter ()}
               unsigned char filter,
               int length,
               array< double >^{\wedge} parameters )
11.69.2.13 FeedbackSetLogic() void FeedbackSetLogic (
               unsigned short position,
               array< unsigned short >^{\wedge} sourcechannel,
               unsigned short resultchannel,
               unsigned int lookup )
11.69.2.14 FeedbackSetMkFilter() void FeedbackSetMkFilter (
               unsigned char filter,
               String^ filtertype,
               double cheb_ribble,
               String^{\wedge} passtype,
               int order,
               double alpha1,
               double alpha2 )
\textbf{11.69.2.15} \quad \textbf{FeedbackSetNumberOfLogics()} \quad \texttt{void FeedbackSetNumberOfLogics} \quad \textbf{(}
               unsigned short number )
11.69.2.16 FeedbackSetNumberOfRateCounter() void FeedbackSetNumberOfRateCounter (
               unsigned short number )
```

```
11.69.2.17 FeedbackSetNumberOfRateDetectors() void FeedbackSetNumberOfRateDetectors (
             unsigned short number )
11.69.2.18 FeedbackSetNumberOfSpikeDetectors() void FeedbackSetNumberOfSpikeDetectors (
             unsigned short number )
11.69.2.19 FeedbackSetNumberOfTriggers() void FeedbackSetNumberOfTriggers (
             unsigned short number )
11.69.2.20 FeedbackSetRateCounter() void FeedbackSetRateCounter (
             unsigned short position,
             unsigned short sourcechannel,
             unsigned short resultchannel )
11.69.2.21 FeedbackSetRateDetector() void FeedbackSetRateDetector (
             unsigned short position,
             unsigned short resultchannel,
             unsigned short trigger,
             unsigned short totzeit,
             unsigned short pulses,
             unsigned int duration1,
             unsigned int duration2 )
11.69.2.22 FeedbackSetSpikeDetectorThreshold() void FeedbackSetSpikeDetectorThreshold (
             unsigned short position,
             unsigned short sourcechannel,
             unsigned short resultchannel,
             unsigned short trigger,
             unsigned short totzeit,
             int threshold1,
             int threshold2,
             short slope )
11.69.2.23 FeedbackSetTrigger() void FeedbackSetTrigger (
             unsigned short position,
             unsigned short sourcechannel,
             unsigned short resultchannel,
             unsigned short trigger,
             unsigned short totzeit )
```

```
11.69.2.24 FeedbackSetTriggerTotzeitFactor() void FeedbackSetTriggerTotzeitFactor (
    unsigned short trigger_factor,
    unsigned short totzeit_factor)
```

11.70 CMealmpedanceDeviceNet Class Reference

Inheritance diagram for CMealmpedanceDeviceNet:

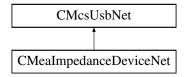

Public Member Functions

- CMealmpedanceDeviceNet ()
- ∼CMeaImpedanceDeviceNet ()
- virtual void StartMeasurement (unsigned short channel)
- virtual unsigned short GetReady ()
- virtual unsigned short GetArraySize ()
- virtual array< unsigned short > ^ GetResult ()
- unsigned short GetAdapterCode ()
- virtual unsigned int GetImpedanceTestFrequency ()
- virtual void SetImpedanceTestFrequency (unsigned int TestFrequency_Hertz)

Additional Inherited Members

11.70.1 Constructor & Destructor Documentation

```
11.70.1.1 CMealmpedanceDeviceNet() CMealmpedanceDeviceNet ( )
```

11.70.1.2 ~ CMealmpedanceDeviceNet() ~ CMealmpedanceDeviceNet ()


11.70.2 Member Function Documentation

 ${\bf 11.70.2.1} \quad {\bf GetAdapterCode()} \quad {\tt unsigned short GetAdapterCode ()}$

11.71 CMeasureTableDeviceNet Class Reference

CMeasureTableDeviceNet is the to control the MCS HLA device

Inheritance diagram for CMeasureTableDeviceNet:

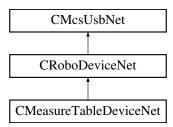

Public Member Functions

CMeasureTableDeviceNet (void)

Properties

• CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.71.1 Detailed Description

CMeasureTableDeviceNet is the to control the MCS HLA device

11.71.2 Constructor & Destructor Documentation

```
11.71.2.1 CMeasureTableDeviceNet() CMeasureTableDeviceNet ( void )
```


11.71.3 Property Documentation

```
11.71.3.1 Sensor CMcsBus_SensorNet^ Sensor [get]
```

11.72 CMeaSwitchDeviceNet Class Reference

The class to control the USB-MEA-Switch.

Inheritance diagram for CMeaSwitchDeviceNet:

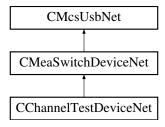

Public Member Functions

CMeaSwitchDeviceNet ()

Constructor.

∼CMeaSwitchDeviceNet ()

Destructor.

• unsigned short GetNumber ()

Gets the number of boards in the device.

array< unsigned char > ^ GetPattern ()

Gets the pattern of the switches that are currently set in the device as char array.

array< bool > ^ GetPatternBool ()

Gets the pattern of the switches that are currently set in he device as bools.

void SetPattern (array< unsigned char >^ pattern)

Sets the pattern of switches from a char array.

void SetPatternBool (array< bool >^ pattern)

Sets the pattern of switches from a.

Additional Inherited Members

11.72.1 Detailed Description

The class to control the USB-MEA-Switch.

This class controls the settings of the USB-MEA-Switch. The box has two inputs for signals from a MEA amplifier. Each of the 64 outputs can be connected to one of the MEAs at the same channel.

11.72.2 Constructor & Destructor Documentation

11.72.2.1 CMeaSwitchDeviceNet() CMeaSwitchDeviceNet ()

Constructor.

11.72.2.2 ~ CMeaSwitchDeviceNet() ~ CMeaSwitchDeviceNet ()

Destructor.

11.72.3 Member Function Documentation

11.72.3.1 GetNumber() unsigned short GetNumber ()

Gets the number of boards in the device.

The MEA-Switch are delivered with 64 or 128 channels

11.72.3.2 GetPattern() array<unsigned char> $^{\land}$ GetPattern ()

Gets the pattern of the switches that are currently set in the device as char array.

11.72.3.3 GetPatternBool() array<bool> ^ GetPatternBool ()

Gets the pattern of the switches that are currently set in he device as bools.

```
11.72.3.4 SetPattern() void SetPattern ( array < unsigned char >^{\wedge} pattern )
```

Sets the pattern of switches from a char array.


```
11.72.3.5 SetPatternBool() void SetPatternBool ( array < bool >^n pattern )
```

Sets the pattern of switches from a.

11.73 CMeaUSBDeviceNet Class Reference

Class for data acquisition via ME and MEA USB amplifiers

Inheritance diagram for CMeaUSBDeviceNet:

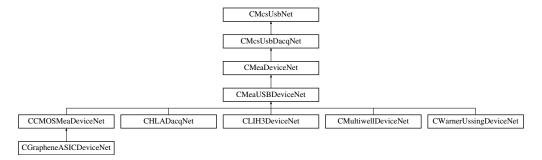

Public Member Functions

- CMeaUSBDeviceNet (OnChannelData^ channelData, OnError^ error)

 Initializes a new instance of CMeaDeviceNet class.
- CMeaUSBDeviceNet ()

Initializes a new instance of CMeaDeviceNet class.

- ∼CMeaUSBDeviceNet ()
- void SetVoltageRangeForIndex (uint16_t index, uint32_t range)

Additional Inherited Members

11.73.1 Detailed Description

Class for data acquisition via ME and MEA USB amplifiers

11.73.2 Constructor & Destructor Documentation

```
11.73.2.1 CMeaUSBDeviceNet() [1/2] CMeaUSBDeviceNet (
OnChannelData^ channelData,
OnError^ error)
```

Initializes a new instance of CMeaDeviceNet class.

chambolbata Handior to dan whom how data is available.	channelData	Handler to call when new data is available.
--	-------------	---

Parameters

```
error Handler to call when an error occurs.
```

11.73.2.2 CMeaUSBDeviceNet() [2/2] CMeaUSBDeviceNet ()


Initializes a new instance of CMeaDeviceNet class.

11.73.2.3 ~ CMeaUSBDeviceNet() ~ CMeaUSBDeviceNet ()

11.73.3 Member Function Documentation

11.74 CMeFunctionNet Class Reference

Inheritance diagram for CMeFunctionNet:

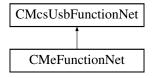

Public Member Functions

CMeFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] meFunctionPointer←
 Container)

Initializes a new instance of the CDacCalibrationFunctionNet class.

- CMeFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CMeFunctionNet (void)
- !CMeFunctionNet (void)
- void SetTransformer (unsigned int index, bool onoff)

Additional Inherited Members

11.74.1 Detailed Description

11.74.2 Constructor & Destructor Documentation

```
11.74.2.1 CMeFunctionNet() [1/2] CMeFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ meFunctionPointerContainer)
```

Initializes a new instance of the CDacCalibrationFunctionNet class.

```
11.74.2.2 CMeFunctionNet() [2/2] CMeFunctionNet ( CMcsUsbNet^{\land} mcsusb )
```


```
11.74.2.3 ~CMeFunctionNet() virtual ~CMeFunctionNet ( void ) [virtual]
```

11.74.3 Member Function Documentation

11.75 CMultiBatteryChargerDeviceNet Class Reference

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

Inheritance diagram for CMultiBatteryChargerDeviceNet:

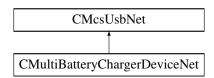

Public Member Functions

CMultiBatteryChargerDeviceNet ()

Initializes a new instance of the CMultiBatteryChargerDeviceNet class.

- virtual \sim CMultiBatteryChargerDeviceNet ()
- !CMultiBatteryChargerDeviceNet ()
- uint32_t GetChargeCurrent (uint32_t NrChannel)

gets the charge current; unit: mA

uint32_t GetDischargeCurrent (uint32_t NrChannel)

gets the discharge current; unit: mA

• void SetDischargeCurrentSetPoint (uint32_t NrChannel, uint32_t DischargeCurrent_mA)

sets the setpoint for the discharge current; unit: mA

uint32 t GetDischargeCurrentSetPoint (uint32 t NrChannel)

gets the setpoint for the discharge current; unit: mA

void SetFinalDischargeVoltage (uint32_t NrChannel, uint32_t FinalDischargeVoltage_mV)

sets the final discharge voltage; unit: mV

uint32_t GetFinalDischargeVoltage (uint32_t NrChannel)

gets the final discharge voltage; unit: mV

uint32_t GetDischargeCapacity (uint32_t NrChannel)

gets the discharge capacity; unit: μAh

uint32_t GetChargeCapacity (uint32_t NrChannel)

gets the charge capacity; unit: μAh

• uint32 t GetBatteryVoltage (uint32 t NrChannel)

gets the battery voltage; unit: mV

uint32_t GetChannels ()

gets number of channels

void SetRatedCapacityVolatile (uint32 t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)

sets the rated capacity (i.e. charge current) without storing it persistently

void SetChargingMode (uint32_t NrChannel, MbcChargingModeEnumNet NewOperatingMode)

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

MbcChargingModeEnumNet GetChargingMode (uint32 t NrChannel)

gets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

• MbcChannelStateEnumNet GetChannelState (uint32_t NrChannel)

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

void CapacityTest (uint32_t NrChannel)

start capacity test on channel

void ChannelReset (uint32 t NrChannel)

cancel charging and capacity test functions; check if battery is connected

void SetChargingPCoefficient (uint32_t pCoefficient)

sets the p-coefficient for charging in mA/V / nominal charging current

uint32 t GetChargingPCoefficient ()

gets the p-coefficient for charging in mA/V / nominal charging current

void SetRatedCapacity (uint32_t NrChannel, MbcRatedCapacityEnumNet NewRatedCapacity)

sets the rated capacity

MbcRatedCapacityEnumNet GetRatedCapacity (uint32 t NrChannel)

gets the rated capacity

Additional Inherited Members

11.75.1 Detailed Description

CMultiBatteryChargerDeviceNet is the class to access the MBC-08 device.

11.75.2 Constructor & Destructor Documentation

```
11.75.2.1 CMultiBatteryChargerDeviceNet() CMultiBatteryChargerDeviceNet ( )
```

Initializes a new instance of the CMultiBatteryChargerDeviceNet class.

```
11.75.2.2 ~CMultiBatteryChargerDeviceNet() virtual ~CMultiBatteryChargerDeviceNet ( ) [virtual]
```

```
11.75.2.3 "!CMultiBatteryChargerDeviceNet() !CMultiBatteryChargerDeviceNet ( )
```

11.75.3 Member Function Documentation

start capacity test on channel

Parameters

NrChannel the channel number

```
11.75.3.2 ChannelReset() void ChannelReset ( uint32_t NrChannel)
```

cancel charging and capacity test functions; check if battery is connected

Parameters

NrChannel the channel number

11.75.3.3 GetBatteryVoltage() uint32_t GetBatteryVoltage (uint32_t NrChannel)

gets the battery voltage; unit: mV

Parameters

NrChannel the channel number

Returns

the battery voltage in mV

11.75.3.4 GetChannels() uint32_t GetChannels ()

gets number of channels

Returns

number of channels

11.75.3.5 **GetChannelState()** MbcChannelStateEnumNet GetChannelState (uint32_t NrChannel)

gets the channel state: IdleNoBattery, IdleChargeFinished, CapacityTestPreCharge, CapacityTestDischarge, StorageCharge, LowCurrentCharge, HighCurrentCharge

Parameters

NrChannel	the channel number

Returns

the current state

11.75.3.6 GetChargeCapacity() uint32_t GetChargeCapacity (uint32_t NrChannel)

gets the charge capacity; unit: µAh

NrChannel the channel number

Returns

the capacity in uAh

gets the charge current; unit: mA

Parameters

NrChannel the channel number

Returns

the measured charge current in mA

11.75.3.8 GetChargingMode() MbcChargingModeEnumNet GetChargingMode (uint32_t NrChannel)

 $gets\ the\ charging\ mode:\ Storage Charge,\ Low Current Charge\ and\ High Current Charge$

Parameters

NrChannel	the channel number
INICHAHILE	lile channel number

Returns

the charging mode

11.75.3.9 **GetChargingPCoefficient()** uint32_t GetChargingPCoefficient ()

gets the p-coefficient for charging in mA/V / nominal charging current

Returns

the p-coefficient

```
11.75.3.10 GetDischargeCapacity() uint32_t GetDischargeCapacity ( uint32_t NrChannel )
```

gets the discharge capacity; unit: μAh

Parameters

NrChannel the channel number

Returns

the capacity in uAh

11.75.3.11 GetDischargeCurrent() uint32_t GetDischargeCurrent (uint32_t NrChannel)

gets the discharge current; unit: mA

Parameters

NrChannel the channel nu

Returns

the measured discharge current in mA

11.75.3.12 GetDischargeCurrentSetPoint() uint32_t GetDischargeCurrentSetPoint (uint32_t NrChannel)

gets the setpoint for the discharge current; unit: mA

Parameters

NrChannel the channel number

Returns

the discharge current in mA

```
11.75.3.13 GetFinalDischargeVoltage() uint32_t GetFinalDischargeVoltage ( uint32_t NrChannel )
```

gets the final discharge voltage; unit: mV

NrChannel	the channel number
-----------	--------------------

Returns

the battery voltage in mV at the end of discharge

11.75.3.14 GetRatedCapacity() MbcRatedCapacityEnumNet GetRatedCapacity (uint32_t NrChannel)

gets the rated capacity

Parameters

NrChannel	the channel number
-----------	--------------------

Returns

the capacity

sets the charging mode: StorageCharge, LowCurrentCharge and HighCurrentCharge

Parameters

NrChannel	the channel number
NewOperatingMode	the charging mode

11.75.3.16 SetChargingPCoefficient() void SetChargingPCoefficient (uint32_t pCoefficient)

sets the p-coefficient for charging in mA/V / nominal charging current

Parameters

-C#:-:+	the meantining
pCoeπicient	the p-coefficient

sets the setpoint for the discharge current; unit: mA

Parameters

NrChannel	the channel number
DischargeCurrent_mA	the discharge current in mA

sets the final discharge voltage; unit: mV

Parameters

NrChannel	the channel number
FinalDischargeVoltage_mV	the battery voltage in mV at the end of discharge

sets the rated capacity

Parameters

NrChannel	the channel number
NewRatedCapacity	the capacity

sets the rated capacity (i.e. charge current) without storing it persistently


Parameters

NrChannel	the channel number
NewRatedCapacity	the capacity

11.76 CMultiwellCallbackFunctionNet Class Reference

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

Inheritance diagram for CMultiwellCallbackFunctionNet:

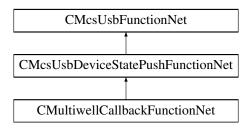

Public Member Functions

- delegate void OnGetPlateClampStateByHeadstage (uint32_t Headstage, PlateClampEnumNet plateState)
- CMultiwellCallbackFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] pMultiwell CallbackFunctionPointerContainer)

Initializes a new instance of the CMultiwellCallbackFunctionNet class.

- CMultiwellCallbackFunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CMultiwellCallbackFunctionNet ()
- !CMultiwellCallbackFunctionNet ()
- PlateClampEnumNet GetPlateClampStateByHeadstage (uint32_t Headstage)

Gets the state of the plate

Events

• OnGetPlateClampStateByHeadstage^ GetPlateClampStateByHeadstageEvent [add, remove, raise]

Event fires when the plate state for the headstage number has changed

Additional Inherited Members

11.76.1 Detailed Description

CMultiwellCallbackFunctionNet is the class to access the Multiwell-Mini-Stimulator

11.76.2 Constructor & Destructor Documentation

```
11.76.2.1 CMultiwellCallbackFunctionNet() [1/2] CMultiwellCallbackFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pMultiwellCallbackFunctionPointerContainer)
```

Initializes a new instance of the CMultiwellCallbackFunctionNet class.

```
11.76.2.2 CMultiwellCallbackFunctionNet() [2/2] CMultiwellCallbackFunctionNet (

CMcsUsbNet^ mcsusb )

11.76.2.3 ~CMultiwellCallbackFunctionNet() virtual ~CMultiwellCallbackFunctionNet ( ) [virtual]

11.76.2.4 "!CMultiwellCallbackFunctionNet() !CMultiwellCallbackFunctionNet ( )

11.76.3 Member Function Documentation
```

11.76.3.1 GetPlateClampStateByHeadstage() PlateClampEnumNet GetPlateClampStateByHeadstage (uint32_t Headstage)

Gets the state of the plate

Parameters

Headstage	The headstage number
-----------	----------------------

Returns

The plate state

11.76.4 Event Documentation


11.76.4.1 GetPlateClampStateByHeadstageEvent OnGetPlateClampStateByHeadstage^ GetPlateClamp↔ StateByHeadstageEvent [add], [remove], [raise]

Event fires when the plate state for the headstage number has changed

11.77 CMultiwellDeviceNet Class Reference

CMultiwellDeviceNet is the class to access the Multiwell device.

Inheritance diagram for CMultiwellDeviceNet:

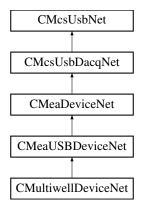

Public Member Functions

CMultiwellDeviceNet ()

Initializes a new instance of the CMultiwellDeviceNet class.

- virtual ∼CMultiwellDeviceNet ()
- !CMultiwellDeviceNet ()
- PlateClampEnumNet GetPlateClampState ()

Gets the state of the Multiwell plate clamp.

• PlateClampEnumNet GetPlateClampState (uint32_t Headstage)

Gets the state of the plate

void OpenPlateClamp ()

Opens the plate clamp.

• void ClosePlateClamp ()

Closes the plate clamp.

void StopPlateClamp ()

Stops the plate clamp movement.

• uint32_t GetPlateClampLockState ()

Gets the state of the plate clamp lock.

• void LockPlateClamp ()

Locks the plate clamp.

void UnlockPlateClamp ()

Unlocks the plate clamp.

MultiwellPlateTypeEnumNet GetPlateType ()

Gets the plate type.

MultiwellPlateTypeEnumNet GetPlateType (uint32_t Headstage)

Gets the plate type.

void SetPlateType (MultiwellPlateTypeEnumNet plateType)

Sets the plate type.

void SetPlateType (uint32_t Headstage, MultiwellPlateTypeEnumNet plateType)

Sets the plate type.

void SetPlateMux (uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

void SetPlateMux (uint32_t Headstage, uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

• uint32 t GetPlateMux ()

Gets the selected quarter of the electrodes on a high density Multiwell plate.

uint32_t GetPlateMux (uint32_t Headstage)

Gets the selected quarter of the electrodes on a high density Multiwell plate.

bool IsPlateTypeValid ()

Checks whether the plate type is valid, meaning all pins have contact.

bool IsPlateTypeValid (uint32_t Headstage)

Checks whether the plate type is valid, meaning all pins have contact.

void SetPowerMuxPlate (uint32_t Headstage, bool powerOn)

On the Multiwell Mini device, turn Power to the MUX Plate On or Off.

bool GetPowerMuxPlate (uint32_t Headstage)

On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.

void SetTouchPadEnable (uint32_t Headstage, bool Enable)

Enables or disables manual opening/closing of plate clamp via touch pad.

bool GetTouchPadEnable (uint32_t Headstage)

Manual opening/closing of plate clamp via touch pad can be disabled.

void SetVolatileClampOffset (uint32_t Headstage, int32_t CoverLipThickness_um)

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

• int32_t GetVolatileClampOffset (uint32_t Headstage)

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Additional Inherited Members

11.77.1 Detailed Description

CMultiwellDeviceNet is the class to access the Multiwell device.

11.77.2 Constructor & Destructor Documentation

11.77.2.1 CMultiwellDeviceNet() CMultiwellDeviceNet ()

Initializes a new instance of the CMultiwellDeviceNet class.

11.77.2.2 ~ CMultiwellDeviceNet() virtual ~ CMultiwellDeviceNet () [virtual]

11.77.2.3 "!CMultiwellDeviceNet() !CMultiwellDeviceNet ()

11.77.3 Member Function Documentation

11.77.3.1 ClosePlateClamp() void ClosePlateClamp ()

Closes the plate clamp.

11.77.3.2 GetPlateClampLockState() uint32_t GetPlateClampLockState ()

Gets the state of the plate clamp lock.

Returns

the state of the plate lock (unlocked/locked)

11.77.3.3 GetPlateClampState() [1/2] PlateClampEnumNet GetPlateClampState ()

Gets the state of the Multiwell plate clamp.

Returns

the state of the plate clamp (open/closed)

11.77.3.4 GetPlateClampState() [2/2] PlateClampEnumNet GetPlateClampState (uint32_t Headstage)

Gets the state of the plate

Parameters

Headstage	The headstage number
-----------	----------------------

Returns

The plate state

11.77.3.5 GetPlateMux() [1/2] uint32_t GetPlateMux ()

Gets the selected quarter of the electrodes on a high density Multiwell plate.

the selected quarter

```
11.77.3.6 GetPlateMux() [2/2] uint32_t GetPlateMux ( uint32_t Headstage )
```

Gets the selected quarter of the electrodes on a high density Multiwell plate.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

the selected quarter

11.77.3.7 GetPlateType() [1/2] MultiwellPlateTypeEnumNet GetPlateType ()

Gets the plate type.

Returns

the plate type

11.77.3.8 GetPlateType() [2/2] MultiwellPlateTypeEnumNet GetPlateType (uint32_t Headstage)

Gets the plate type.

Parameters

Headstage The he	eadstage to query.
------------------	--------------------

Returns

the plate type

11.77.3.9 **GetPowerMuxPlate()** bool GetPowerMuxPlate (uint32_t *Headstage*)

On the Multiwell Mini device, Query if Power to the MUX Plate is On or Off.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

"true" Power is On, "false" Power is Off

11.77.3.10 GetTouchPadEnable() bool GetTouchPadEnable (uint32_t Headstage)

Manual opening/closing of plate clamp via touch pad can be disabled.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

"true" when plate clamp can be driven manually, otherwise "false".

11.77.3.11 GetVolatileClampOffset() int32_t GetVolatileClampOffset (uint32_t Headstage)

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Parameters

Headstage	The headstage to be affected.

Returns

Clamp distance reduction in um.

11.77.3.12 IsPlateTypeValid() [1/2] bool IsPlateTypeValid ()

Checks whether the plate type is valid, meaning all pins have contact.

Returns

"true" when all pins have contact, otherwise "false".

11.77.3.13 **IsPlateTypeValid()** [2/2] bool IsPlateTypeValid (uint32_t *Headstage*)

Checks whether the plate type is valid, meaning all pins have contact.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

"true" when all pins have contact, otherwise "false".

11.77.3.14 LockPlateClamp() void LockPlateClamp ()

Locks the plate clamp.

11.77.3.15 OpenPlateClamp() void OpenPlateClamp ()

Opens the plate clamp.

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

Headstage	The headstage to query.
muxSelection	the selected quarter

11.77.3.17 SetPlateMux() [2/2] void SetPlateMux (uint32_t muxSelection)

Selects a one quarter of the electrodes on a high density Multiwell plate.

Parameters

muxSelection	the selected quarter

Sets the plate type.

Parameters

```
plateType the plate type
```

Sets the plate type.

Parameters

Headstage	The headstage to query.
plateType	the plate type

On the Multiwell Mini device, turn Power to the MUX Plate On or Off.

Parameters

Headstage	The headstage to query.
powerOn	"true" to turn Power On, "false" to turn Power Off

Enables or disables manual opening/closing of plate clamp via touch pad.

Parameters

Headstage	The headstage to be affected.
Enable	"true" when plate clamp shall be driven manually, otherwise "false".

The distance travelled to clamp the plate can be reduced e.g. to compensate for a cover.

Parameters

Headstage	The headstage to be affected.
CoverLipThickness_um	Clamp distance reduction in um. Range: 200um400um; Typical value: 300um

11.77.3.23 StopPlateClamp() void StopPlateClamp ()

Stops the plate clamp movement.


11.77.3.24 UnlockPlateClamp() void UnlockPlateClamp ()

Unlocks the plate clamp.

11.78 CMultiwellOptoStimFunctionNet Class Reference

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

Inheritance diagram for CMultiwellOptoStimFunctionNet:

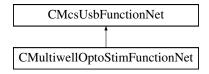

Public Member Functions

CMultiwellOptoStimFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pMultiwell← OptoStimFunctionPointerContainer)

Initializes a new instance of the CMultiwellOptoStimFunctionNet class.

- CMultiwellOptoStimFunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CMultiwellOptoStimFunctionNet ()
- !CMultiwellOptoStimFunctionNet ()
- uint32_t GetWaveLengthInNanometer (uint16_t channel)
- uint32 t GetAbsMaxCurrentInMicroAmp (uint16 t channel)
- uint32_t GetMaxDurationHighCurrentInMicroSec (uint16_t channel)

- uint32_t GetMaxDutyCycleHighCurrent (uint16_t channel)
- uint32 t GetPermanentCurrentInMicroAmp (uint16 t channel)
- uint32_t GetColorRgb (uint16_t channel)
- String \(^\) GetColorStr (uint16 t channel)
- void SetWaveLengthInNanometer (uint16 t channel, uint32 t WaveLength nm)
- void SetAbsMaxCurrentInMicroAmp (uint16_t channel, uint32_t AbsoluteMaxCurrent_uA)
- void SetMaxDurationHighCurrentInMicroSec (uint16_t channel, uint32_t AbsoluteMaxDuration_us)
- void SetMaxDutyCycleHighCurrent (uint16_t channel, uint32_t MaxDutyCycleHighCurrent)
- void SetPermanentCurrentInMicroAmp (uint16_t channel, uint32_t PermanentCurrent_uA)
- void SetColorRgb (uint16_t channel, uint32_t ColorRGB)
- void SetColorStr (uint16 t channel, String[^] ColorString)

Additional Inherited Members

11.78.1 Detailed Description

CMultiwellOptoStimFunctionNet is the class to access the optical properties of the Multiwell Optostim device

11.78.2 Constructor & Destructor Documentation

```
11.78.2.1 CMultiwellOptoStimFunctionNet() [1/2] CMultiwellOptoStimFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pMultiwellOptoStimFunctionPointerContainer )
```

Initializes a new instance of the CMultiwellOptoStimFunctionNet class.

```
11.78.2.2 CMultiwellOptoStimFunctionNet() [2/2] CMultiwellOptoStimFunctionNet (
CMcsUsbNet^ mcsusb)
```

```
11.78.2.3 ~CMultiwellOptoStimFunctionNet() virtual ~CMultiwellOptoStimFunctionNet ( ) [virtual]
```

```
11.78.2.4 "!CMultiwellOptoStimFunctionNet() !CMultiwellOptoStimFunctionNet ( )
```

11.78.3 Member Function Documentation

```
11.78.3.1 GetAbsMaxCurrentInMicroAmp() uint32_t GetAbsMaxCurrentInMicroAmp ( uint16_t channel)
```

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

absolute max. current; unit: uA

11.78.3.2 GetColorRgb() uint32_t GetColorRgb (uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

RGB-value of LED color

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

LED color as string

11.78.3.4 GetMaxDurationHighCurrentInMicroSec() uint32_t GetMaxDurationHighCurrentInMicroSec (uint16_t channel)

Parameters

channel the (analog) channel i	number
--------------------------------	--------

Returns

max. duration the LED can stand the abs. max current; unit: us

11.78.3.5 GetMaxDutyCycleHighCurrent() uint32_t GetMaxDutyCycleHighCurrent (

uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

max. duty cycle at max. current; unit: 100*%

11.78.3.6 GetPermanentCurrentInMicroAmp() uint32_t GetPermanentCurrentInMicroAmp (

uint16_t channel)

Parameters

channel	the (analog) channel number
---------	-----------------------------

Returns

max. current the LED can stand when always switched on; unit: uA

11.78.3.7 GetWaveLengthInNanometer() uint32_t GetWaveLengthInNanometer (

uint16_t channel)

Parameters

channel the (analog) channel numb	er
-----------------------------------	----

Returns

wavelength of this channel's LEDs; unit: nm

11.78.3.8 SetAbsMaxCurrentInMicroAmp() void SetAbsMaxCurrentInMicroAmp (

```
uint16_t channel,
uint32_t AbsoluteMaxCurrent_uA )
```

Parameters

channel	the (analog) channel number
AbsoluteMaxCurrent_uA	absolute max. current; unit: uA

Parameters

channel	the (analog) channel number
ColorRGB	RGB-value of LED color

Parameters

channel	the (analog) channel number
ColorString	LED color as string

Parameters

channel	the (analog) channel number
AbsoluteMaxDuration_us	max. duration the LED can stand the abs. max current; unit: us

Parameters

channel	the (analog) channel number
MaxDutyCycleHighCurrent	max. duty cycle at max. current; unit: 100*%

11.78.3.13 SetPermanentCurrentInMicroAmp() void SetPermanentCurrentInMicroAmp (

```
uint16_t channel,
uint32_t PermanentCurrent_uA )
```

Parameters


channel	the (analog) channel number	
PermanentCurrent_uA	max. current the LED can stand when always switched on; unit: uA	

Parameters

channel	the (analog) channel number
WaveLength_nm	wavelength of this channel's LEDs; unit: nm

11.79 CNF_GenDeviceNet Class Reference

Inheritance diagram for CNF_GenDeviceNet:

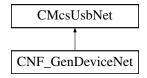

Public Member Functions

- CNF_GenDeviceNet (void)
- \sim CNF_GenDeviceNet (void)
- void Set_Values (unsigned int frequency, unsigned int amplitude)

Additional Inherited Members

11.79.1 Constructor & Destructor Documentation

11.79.1.1 CNF_GenDeviceNet() CNF_GenDeviceNet (void)


```
11.79.1.2 ~CNF_GenDeviceNet() ~CNF_GenDeviceNet (
```

11.79.2 Member Function Documentation

```
11.79.2.1 Set_Values() void Set_Values (
          unsigned int frequency,
          unsigned int amplitude )
```

11.80 COctoPotDeviceNet Class Reference

Inheritance diagram for COctoPotDeviceNet:

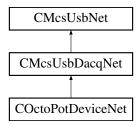

Public Member Functions

- COctoPotDeviceNet (void)
- COctoPotDeviceNet (OnChannelData^ channelData, OnError^ error)
- uint32_t SetOutputRate (uint32_t rate)
- uint32_t SetBathclamp (unsigned int block, bool enable)
- uint32_t SetDacValue (int channel, int value)
- uint32 t SetDacAutoControl (unsigned int channel)
- uint32_t SetPidParameter (unsigned int channel, int const_p, int const_i, int shift_p, int shift_i)
- uint32_t SetRampParameter (unsigned int channel, int start, int min, int max, int slope, int slope2, int pause, unsigned int samples)
- uint32_t RampStart (int channelmap)
- uint32 t SetSineParameter (unsigned int channel, int amplitude)
- uint32_t SineStart (int channelmap)
- · uint32_t SetPatternListEntry (unsigned int channel, unsigned int position, unsigned int duration, int value)
- uint32_t PatternListStart (int channelmap)
- uint32_t SetAdcOffset (unsigned int channel, int offset)
- uint32 t SetDacOffset (unsigned int channel, int offset)
- uint32_t ResetAdcOffset (unsigned int channel)
- uint32_t ResetDacOffset (unsigned int channel)
- uint32 t BurnAdcOffset ()
- uint32_t BurnDacOffset ()
- uint32_t GetAdcOffset (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- uint32_t GetDacOffset (unsigned int channel, [System::Runtime::InteropServices::Out] int ^ offset)
- uint32_t SetAmplificationSwitch (unsigned int channel, unsigned int state)
- uint32_t SetChannelSwitch (unsigned int channel, unsigned int state)
- uint32_t SetNumberOfChannels (unsigned int NumberOfChannels)
- uint32 t EnableDigitalIn (bool enable)
- uint32_t EnableTimestamp (bool enable)
- uint32_t EnableChecksum (bool enable)

Additional Inherited Members

11.80.1 Constructor & Destructor Documentation

```
11.80.1.1 COctoPotDeviceNet() [1/2] COctoPotDeviceNet (
             void )
11.80.1.2 COctoPotDeviceNet() [2/2] COctoPotDeviceNet (
             OnChannelData^ channelData,
             OnError<sup>∧</sup> error )
11.80.2 Member Function Documentation
11.80.2.1 BurnAdcOffset() uint32_t BurnAdcOffset ( )
11.80.2.2 BurnDacOffset() uint32_t BurnDacOffset ()
11.80.2.3 EnableChecksum() uint32_t EnableChecksum (
             bool enable )
11.80.2.4 EnableDigitalIn() uint32_t EnableDigitalIn (
             bool enable )
11.80.2.5 EnableTimestamp() uint32_t EnableTimestamp (
             bool enable )
11.80.2.6 GetAdcOffset() uint32_t GetAdcOffset (
             unsigned int channel,
             [System::Runtime::InteropServices::Out] int ^{\wedge} offset )
```


```
11.80.2.7 GetDacOffset() uint32_t GetDacOffset (
             unsigned int channel,
             [System::Runtime::InteropServices::Out] int ^{\land} offset )
11.80.2.8 PatternListStart() uint32_t PatternListStart (
             int channelmap )
11.80.2.9 RampStart() uint32_t RampStart (
             int channelmap )
11.80.2.10 ResetAdcOffset() uint32_t ResetAdcOffset (
             unsigned int channel )
11.80.2.11 ResetDacOffset() uint32_t ResetDacOffset (
             unsigned int channel )
11.80.2.12 SetAdcOffset() uint32_t SetAdcOffset (
             unsigned int channel,
             int offset )
11.80.2.13 SetAmplificationSwitch() uint32_t SetAmplificationSwitch (
             unsigned int channel,
             unsigned int state )
11.80.2.14 SetBathclamp() uint32_t SetBathclamp (
             unsigned int block,
             bool enable )
11.80.2.15 SetChannelSwitch() uint32_t SetChannelSwitch (
             unsigned int channel,
             unsigned int state )
```

```
11.80.2.16 SetDacAutoControl() uint32_t SetDacAutoControl (
             unsigned int channel )
11.80.2.17 SetDacOffset() uint32_t SetDacOffset (
             unsigned int channel,
             int offset )
11.80.2.18 SetDacValue() uint32_t SetDacValue (
             int channel,
             int value )
11.80.2.19 SetNumberOfChannels() uint32_t SetNumberOfChannels (
             unsigned int NumberOfChannels )
11.80.2.20 SetOutputRate() uint32_t SetOutputRate (
             uint32_t rate )
11.80.2.21 SetPatternListEntry() uint32_t SetPatternListEntry (
             unsigned int channel,
             unsigned int position,
             unsigned int duration,
             int value )
11.80.2.22 SetPidParameter() uint32_t SetPidParameter (
             unsigned int channel,
             int const_p,
             int const_i,
             int shift_p,
             int shift_i )
11.80.2.23 SetRampParameter() uint32_t SetRampParameter (
             unsigned int channel,
             int start,
             int min,
             int max,
             int slope,
             int slope2,
             int pause,
             unsigned int samples )
```

```
11.80.2.24 SetSineParameter() uint32_t SetSineParameter (
          unsigned int channel,
          int amplitude )
11.80.2.25 SineStart() uint32_t SineStart (
          int channelmap )
```

11.81 COkuvisionStimulatorDeviceNet Class Reference

Inheritance diagram for COkuvisionStimulatorDeviceNet:

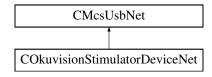

Public Member Functions

- · COkuvisionStimulatorDeviceNet (void)
- ~COkuvisionStimulatorDeviceNet (void)
- void SetPulseform (int channel, int current, int pulsewidth, int periode, int duration)
- void GetPulseform (int channel, [System::Runtime::InteropServices::Out] int% current, [System::Runtime
 ::InteropServices::Out] int% pulsewidth, [System::Runtime::InteropServices::Out] int% periode, [System::
 Runtime::InteropServices::Out] int% duration)
- void SetMaxPower (int channel, int power)
- int GetMaxPower (int channel)
- void SetMaxVoltage (int channel, int voltage)
- int GetMaxVoltage (int channel)
- void SetCheckVoltage (int channel, int voltage)
- int GetCheckVoltage (int channel)
- int GetVoltage (int channel)
- void SetDACOffset (int channel, int part, int offset)
- · int GetDACOffset (int channel, int part)
- void SetRTC (uint8_t year, uint8_t month, uint8_t day, uint8_t hour, uint8_t minute, uint8_t second)
- void GetRTC ([System::Runtime::InteropServices::Out] uint8_t% year, [System::Runtime::InteropServices
 ::Out] uint8_t% month, [System::Runtime::InteropServices::Out] uint8_t% day, [System::Runtime::Interop
 Services::Out] uint8_t% hour, [System::Runtime::InteropServices::Out] uint8_t% minute, [System::Runtime
 ::InteropServices::Out] uint8_t% second)
- void SetRTC (DateTime timestamp)
- DateTime GetRTC ()
- void GetStimulatorStatus ([System::Runtime::InteropServices::Out] int% startstop, [System::Runtime::
 — InteropServices::Out] int% last_error, [System::Runtime::InteropServices::Out] int% battery_status)
- void SetCurrentFactor (int channel, int factor)
- int GetCurrentFactor (int channel)

Additional Inherited Members

11.81.1 Constructor & Destructor Documentation

```
11.81.1.1 COkuvisionStimulatorDeviceNet() COkuvisionStimulatorDeviceNet (
            void )
11.81.1.2 ~COkuvisionStimulatorDeviceNet() ~COkuvisionStimulatorDeviceNet (
             void )
11.81.2 Member Function Documentation
11.81.2.1 GetCheckVoltage() int GetCheckVoltage (
            int channel )
11.81.2.2 GetCurrentFactor() int GetCurrentFactor (
            int channel )
11.81.2.3 GetDACOffset() int GetDACOffset (
            int channel,
            int part )
11.81.2.4 GetMaxPower() int GetMaxPower (
            int channel )
11.81.2.5 GetMaxVoltage() int GetMaxVoltage (
            int channel )
```


```
11.81.2.6 GetPulseform() void GetPulseform (
             int channel,
             [System::Runtime::InteropServices::Out] int% current,
             [{\tt System::Runtime::InteropServices::Out}] \  \, {\tt int\$} \  \, {\tt pulsewidth},
             [System::Runtime::InteropServices::Out] int% periode,
             [System::Runtime::InteropServices::Out] int% duration )
11.81.2.7 GetRTC() [1/2] DateTime GetRTC ( )
11.81.2.8 GetRTC() [2/2] void GetRTC (
             [System::Runtime::InteropServices::Out] uint8_t% year,
             [System::Runtime::InteropServices::Out] uint8_t% month,
             [System::Runtime::InteropServices::Out] uint8_t% day,
             [System::Runtime::InteropServices::Out] uint8_t% hour,
             [System::Runtime::InteropServices::Out] uint8_t% minute,
             [System::Runtime::InteropServices::Out] uint8_t% second )
11.81.2.9 GetStimulatorStatus() void GetStimulatorStatus (
             [System::Runtime::InteropServices::Out] int% startstop,
             [System::Runtime::InteropServices::Out] int% last_error,
             [System::Runtime::InteropServices::Out] int% battery_status )
11.81.2.10 GetVoltage() int GetVoltage (
             int channel )
11.81.2.11 SetCheckVoltage() void SetCheckVoltage (
             int channel,
             int voltage )
11.81.2.12 SetCurrentFactor() void SetCurrentFactor (
             int channel,
             int factor )
```

```
11.81.2.13 SetDACOffset() void SetDACOffset (
             int channel,
             int part,
             int offset )
11.81.2.14 SetMaxPower() void SetMaxPower (
             int channel,
             int power )
11.81.2.15 SetMaxVoltage() void SetMaxVoltage (
             int channel,
             int voltage )
11.81.2.16 SetPulseform() void SetPulseform (
             int channel,
             int current,
             int pulsewidth,
             int periode,
             int duration )
11.81.2.17 SetRTC() [1/2] void SetRTC (
             \texttt{DateTime} \ \textit{timestamp} \ )
11.81.2.18 SetRTC() [2/2] void SetRTC (
             uint8_t year,
             uint8_t month,
             uint8_t day,
             uint8_t hour,
             uint8_t minute,
             uint8_t second )
```

11.82 CPatchServerDeviceNet Class Reference

CPatchServerDeviceNet is the class to control the MCS PatchServer device

Inheritance diagram for CPatchServerDeviceNet:

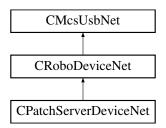

Public Member Functions

CPatchServerDeviceNet (void)

Properties

• CMcsBus_SensorNet^ Sensor [get]

Additional Inherited Members

11.82.1 Detailed Description

CPatchServerDeviceNet is the class to control the MCS PatchServer device

11.82.2 Constructor & Destructor Documentation


```
11.82.2.1 CPatchServerDeviceNet() CPatchServerDeviceNet (
```

11.82.3 Property Documentation

```
11.82.3.1 Sensor CMcsBus_SensorNet^ Sensor [get]
```

11.83 CPathIdentDeviceNet Class Reference

Inheritance diagram for CPathIdentDeviceNet:

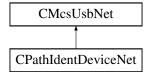

Public Member Functions

- CPathIdentDeviceNet (void)
- ∼CPathIdentDeviceNet (void)
- void Set_Values (unsigned int frequency, unsigned int amplitude)


Additional Inherited Members

11.83.1 Constructor & Destructor Documentation

11.84 CPedoterDeviceNet Class Reference

unsigned int amplitude)

Inheritance diagram for CPedoterDeviceNet:

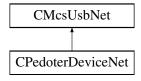

Public Member Functions

• CPedoterDeviceNet ()

Initializes a new instance of the CPedoterDeviceNet class.

- virtual ∼CPedoterDeviceNet ()
- !CPedoterDeviceNet ()
- uint32_t GetCommand (uint16_t Argument)

Get value from the pedoter device

void SetCommand (uint16_t Argument, uint32_t pData)

Set value on the pedoter device

Additional Inherited Members

11.84.1 Detailed Description

11.84.2 Constructor & Destructor Documentation

```
11.84.2.1 CPedoterDeviceNet() CPedoterDeviceNet ()
```

Initializes a new instance of the CPedoterDeviceNet class.

```
11.84.2.2 ~ CPedoterDeviceNet() virtual ~ CPedoterDeviceNet () [virtual]
```

```
11.84.2.3 "!CPedoterDeviceNet() !CPedoterDeviceNet ( )
```

11.84.3 Member Function Documentation

```
11.84.3.1 GetCommand() uint32_t GetCommand ( uint16_t Argument )
```

Get value from the pedoter device

Parameters

```
Argument argument
```

Returns

value

Set value on the pedoter device


Parameters

Argument	argument
pData	value

11.85 CPeristalticPumpDeviceNet Class Reference

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

Inheritance diagram for CPeristalticPumpDeviceNet:

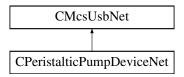

Public Member Functions

- CPeristalticPumpDeviceNet (void)

 Initialize a new instance of the CPeristalticPumpDeviceNet class.
- ~CPeristalticPumpDeviceNet (void)

Properties

• CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

Additional Inherited Members

11.85.1 Detailed Description

CPeristalticPumpDeviceNet is the class to control a Persistaltic Pump.

11.85.2 Constructor & Destructor Documentation

Initialize a new instance of the CPeristalticPumpDeviceNet class.


```
11.85.2.2 ~CPeristalticPumpDeviceNet() ~CPeristalticPumpDeviceNet (
```

11.85.3 Property Documentation

11.85.3.1 McsBus MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.86 CPgaDeviceNet Class Reference

Inheritance diagram for CPgaDeviceNet:

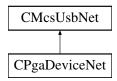

Public Member Functions

- CPgaDeviceNet ()
- ∼CPgaDeviceNet ()
- uint32_t GetNumFrequencyRanges ([System::Runtime::InteropServices::Out]int% numRanges)
- uint32_t GetFrequencyRange (int rangeIndex, [System::Runtime::InteropServices::Out]int% low, [System::

 Runtime::InteropServices::Out]int% high, [System::Runtime::InteropServices::Out]int% channels, [System
 ::Runtime::InteropServices::Out]int% gain)
- uint32_t GetNumAmplifications ([System::Runtime::InteropServices::Out]int% number)
- uint32_t GetAmplification (int index, [System::Runtime::InteropServices::Out]int% amplification, [System::
 Runtime::InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)
- uint32 t DefineNumFrequencyRanges (int rnum)
- uint32 t DefineFrequencyRange (int index, int low, int high, int channels, int gain)
- uint32 t DefineNumAmplifications (int number)
- uint32_t DefineAmplification (int index, int amplification, int poti1, int poti2)
- uint32_t SetGain (int channel, int Gain, int poti1, int poti2)
- uint32_t GetGain (int channel, [System::Runtime::InteropServices::Out]int% Gain, [System::Runtime::← InteropServices::Out]int% poti1, [System::Runtime::InteropServices::Out]int% poti2)
- uint32 t ApplyGains ()

Additional Inherited Members

11.86.1 Constructor & Destructor Documentation

11.86.1.1 CPgaDeviceNet() CPgaDeviceNet ()


```
11.86.1.2 ~CPgaDeviceNet() ~CPgaDeviceNet ()
11.86.2 Member Function Documentation
11.86.2.1 ApplyGains() uint32_t ApplyGains ( )
11.86.2.2 DefineAmplification() uint32_t DefineAmplification (
             int index,
             int amplification,
             int poti1,
             int poti2 )
11.86.2.3 DefineFrequencyRange() uint32_t DefineFrequencyRange (
             int index,
             int low,
             int high,
             int channels,
             int gain )
11.86.2.4 DefineNumAmplifications() uint32_t DefineNumAmplifications (
             int number )
11.86.2.5 DefineNumFrequencyRanges() uint32_t DefineNumFrequencyRanges (
             int rnum )
11.86.2.6 GetAmplification() uint32_t GetAmplification (
             int index,
             [System::Runtime::InteropServices::Out] int% amplification,
             [System::Runtime::InteropServices::Out] int% poti1,
             [System::Runtime::InteropServices::Out] int% poti2 )
```

```
11.86.2.7 GetFrequencyRange() uint32_t GetFrequencyRange (
             int rangeIndex,
             [System::Runtime::InteropServices::Out] int% low,
             [System::Runtime::InteropServices::Out] int% high,
             [System::Runtime::InteropServices::Out] int% channels,
             [System::Runtime::InteropServices::Out] int% gain )
11.86.2.8 GetGain() uint32_t GetGain (
             int channel,
             [System::Runtime::InteropServices::Out] int% Gain,
             [System::Runtime::InteropServices::Out] int% potil,
             [System::Runtime::InteropServices::Out] int% poti2 )
11.86.2.9 GetNumAmplifications() uint32_t GetNumAmplifications (
             [System::Runtime::InteropServices::Out] int% number )
11.86.2.10 GetNumFrequencyRanges() uint32_t GetNumFrequencyRanges (
             [System::Runtime::InteropServices::Out] int% numRanges )
11.86.2.11 SetGain() uint32_t SetGain (
             int channel,
             int Gain,
             int poti1,
             int poti2 )
```

11.87 CPositionIIDeviceNet Class Reference

CPositionIIDeviceNet is the class to control PositionII devices

Inheritance diagram for CPositionIIDeviceNet:

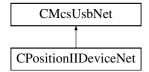

Public Member Functions

CPositionIIDeviceNet ()

Initializes a new instance of the CPositionIIDeviceNet class.

- virtual ∼CPositionIIDeviceNet ()
- !CPositionIIDeviceNet ()
- uint32 t GetCoilCommunication (uint16 t coil)

get if the communication to the coil is working

uint32 t GetOnOff (uint16 t coil)

get if the coil is switched on/off

• void SwitchOnOff (uint16_t coil, uint32_t on)

switched the coild on of

uint32 t GetImplantState (uint16 t coil)

gets the implantat state

uint32_t GetImplantCurrentSetpoint (uint16_t coil)

sets the implant current setpoint

· void SetImplantCurrentSetpoint (uint16 t coil, uint32 t current)

gets the implant current setpoint

• uint32 t GetPowerStrength (uint16 t coil)

sets the power for the trigger pulses

· void SetPowerStrength (uint16_t coil, uint32_t power)

gets the power for the trigger pulses

· uint32 t GetImplantResult (uint16 t coil)

gets the last result of the implant pulse trigger

void GetRTC ([System::Runtime::InteropServices::Out]uint8_t% year, [System::Runtime::InteropServices::Out]uint8_t% month, [System::Runtime::InteropServices::Out]uint8_t% day, [System::Runtime::InteropServices::Out]uint8_t% hour, [System::Runtime::InteropServices::Out]uint8_t% minute, [System::Runtime::InteropServices::Out]uint8_t% second)

Get the RTC

void SetRTC (uint8 t year, uint8 t month, uint8 t day, uint8 t hour, uint8 t minute, uint8 t second)

Set the RTC

uint32_t GetStateDebugData (uint16_t coil)

get the debug queue state

void SetStateDebugData (uint16_t coil, uint32_t state)

clears/starts/stops the debug queue for a certain coil

void GetDebugData (uint16_t coil, [System::Runtime::InteropServices::Out]uint16_t% index, [System::
 Runtime::InteropServices::Out]uint16_t% voltage, [System::Runtime::InteropServices::Out]uint16_t% numberofpulses, [System::Runtime::InteropServices::Out]uint16_t% mediantime)

get the oldest debug entry for a certain coil

• uint32_t GetStateEventData ()

get the event queue state

void SetStateEventData (uint32_t state)

clears/starts/stops the event queue for a certain coil

• void GetEventData ([System::Runtime::InteropServices::Out]uint16_t% index, [System::Runtime::Interop Services::Out]uint8_t% year, [System::Runtime::InteropServices::Out]uint8_t% month, [System::Runtime::InteropServices::Out]uint8_t% day, [System::Runtime::InteropServices::Out]uint8_t% hour, [System::Cout]uint8_t% minute, [System::Runtime::InteropServices::Out]uint8_t% second, [System::Runtime::InteropServices::Out]uint16_t% coil, [System::Runtime::InteropServices::Out]uint16_t% type, [System::Runtime::InteropServices::Out]uint16_t% value)

get the oldest event entry

Properties

• CRFFunctionNet^ RFFunction [get]

Additional Inherited Members

11.87.1 Detailed Description

CPositionIIDeviceNet is the class to control PositionII devices

11.87.2 Constructor & Destructor Documentation

```
11.87.2.1 CPositionIIDeviceNet() CPositionIIDeviceNet ()
```

Initializes a new instance of the CPositionIIDeviceNet class.

```
11.87.2.2 ~CPositionIIDeviceNet() virtual ~CPositionIIDeviceNet ( ) [virtual]
```

```
11.87.2.3 "!CPositionIIDeviceNet() !CPositionIIDeviceNet ( )
```

11.87.3 Member Function Documentation

```
11.87.3.1 GetCoilCommunication() uint32_t GetCoilCommunication ( uint16_t coil)
```

get if the communication to the coil is working

Parameters

coil the coil

Returns

is communicating

11.87.3.2 GetDebugData() void GetDebugData (

```
uint16_t coil,
[System::Runtime::InteropServices::Out] uint16_t% index,
[System::Runtime::InteropServices::Out] uint16_t% voltage,
[System::Runtime::InteropServices::Out] uint16_t% numberofpulses,
[System::Runtime::InteropServices::Out] uint16_t% mediantime)
```

get the oldest debug entry for a certain coil

Parameters

coil	the coil
index	the debug entry index number
voltage	the voltage applied
numberofpulses	the number of pulses detected
mediantime	the median time between pulses

11.87.3.3 GetEventData() void GetEventData (

```
[System::Runtime::InteropServices::Out] uint16_t% index,
[System::Runtime::InteropServices::Out] uint8_t% year,
[System::Runtime::InteropServices::Out] uint8_t% month,
[System::Runtime::InteropServices::Out] uint8_t% day,
[System::Runtime::InteropServices::Out] uint8_t% hour,
[System::Runtime::InteropServices::Out] uint8_t% minute,
[System::Runtime::InteropServices::Out] uint8_t% second,
[System::Runtime::InteropServices::Out] uint16_t% coil,
[System::Runtime::InteropServices::Out] uint16_t% type,
[System::Runtime::InteropServices::Out] uint16_t% value)
```

get the oldest event entry

Parameters

index	the event index number
year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second
coil	the coil
type	the event type
value	the even value

```
11.87.3.4 GetImplantCurrentSetpoint() uint32_t GetImplantCurrentSetpoint ( uint16_t coil )
```

sets the implant current setpoint

Parameters

coil	the coil
------	----------

the current

```
11.87.3.5 GetImplantResult() uint32_t GetImplantResult ( uint16_t coil )
```

gets the last result of the implant pulse trigger

Parameters

```
coil the coil
```

Returns

the result

```
11.87.3.6 GetImplantState() uint32_t GetImplantState ( uint16_t coil )
```

gets the implantat state

Parameters

```
coil the coil
```

Returns

the state

```
11.87.3.7 GetOnOff() uint32_t GetOnOff ( uint16_t coil )
```

get if the coil is switched on/off

Parameters

coil the coil

```
0 = off, 1 = on
```

11.87.3.8 GetPowerStrength() uint32_t GetPowerStrength (uint16_t coil)

sets the power for the trigger pulses

Parameters

coil	the coil
------	----------

Returns

the power

11.87.3.9 GetRTC() void GetRTC (

```
[System::Runtime::InteropServices::Out] uint8_t% year,
[System::Runtime::InteropServices::Out] uint8_t% month,
[System::Runtime::InteropServices::Out] uint8_t% day,
[System::Runtime::InteropServices::Out] uint8_t% hour,
[System::Runtime::InteropServices::Out] uint8_t% minute,
[System::Runtime::InteropServices::Out] uint8_t% second )
```

Get the RTC

Parameters

year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second

11.87.3.10 GetStateDebugData() uint32_t GetStateDebugData (uint16_t coil)

get the debug queue state

Parameters

coil the coil

the state

```
11.87.3.11 GetStateEventData() uint32_t GetStateEventData ( )
```

get the event queue state

Returns

the state

gets the implant current setpoint

Parameters

coil	the coil
current	the current

gets the power for the trigger pulses

Parameters

coil	the coil
power	the power

```
11.87.3.14 SetRTC() void SetRTC (
    uint8_t year,
    uint8_t month,
    uint8_t day,
    uint8_t hour,
    uint8_t minute,
    uint8_t second )
```

Set the RTC

year	the year
month	the month
day	the day
hour	the hour
minute	the minute
second	the second

clears/starts/stops the debug queue for a certain coil

Parameters

coil	the coil
state	clear/start/stop

11.87.3.16 SetStateEventData() void SetStateEventData (uint32_t state)

clears/starts/stops the event queue for a certain coil

Parameters

state	clear/start/stop

```
11.87.3.17 SwitchOnOff() void SwitchOnOff (
    uint16_t coil,
    uint32_t on)
```

switched the coild on of

Parameters

coil	the coil
on	0 = off, 1 = on


11.87.4 Property Documentation

11.87.4.1 RFFunction CRFFunctionNet^ RFFunction [get]

11.88 CPositionImpDeviceNet Class Reference

CPositionImpDeviceNet is the class to access the Position/Imp devices

Inheritance diagram for CPositionImpDeviceNet:

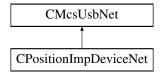

Public Member Functions

CPositionImpDeviceNet ()

Initializes a new instance of the CPositionImpDeviceNet class.

- virtual ∼CPositionImpDeviceNet ()
- !CPositionImpDeviceNet ()
- void ConnectImp (uint32_t id)

Connect to a Imp device with a certain ID

uint32_t ConnectedImp ()

The ID of the connected Imp device

• int32_t GetRFFrequency ()

Gets currently used RF frequency

void SetRFFrequency (int32_t frequency)

Sets the current RF frequency

• uint32_t GetDeviceList (int32_t index)

Gets the device list

• void SetDeviceList (int32_t index, uint32_t id)

Sets the device list

• uint32_t GetImpId ()

Gets the ID of the impedance measure device

void SetImpId (uint32_t id)

Sets the ID of the impedance measure device

Additional Inherited Members

11.88.1 Detailed Description

CPositionImpDeviceNet is the class to access the Position/Imp devices

11.88.2 Constructor & Destructor Documentation

```
11.88.2.1 CPositionImpDeviceNet() CPositionImpDeviceNet ()
```

Initializes a new instance of the CPositionImpDeviceNet class.

```
11.88.2.2 ~CPositionImpDeviceNet() virtual ~CPositionImpDeviceNet () [virtual]
```

```
11.88.2.3 "!CPositionImpDeviceNet() !CPositionImpDeviceNet ( )
```

11.88.3 Member Function Documentation

11.88.3.1 ConnectedImp() uint32_t ConnectedImp ()

The ID of the connected Imp device

Returns

The ID

```
11.88.3.2 ConnectImp() void ConnectImp ( uint32_t id )
```

Connect to a Imp device with a certain ID

Parameters

id The ID

Gets the device list

Returns

the ID

```
11.88.3.4 GetImpId() uint32_t GetImpId ( )
```

Gets the ID of the impedance measure device

Returns

the ID

```
11.88.3.5 GetRFFrequency() int32_t GetRFFrequency ( )
```

Gets currently used RF frequency

Returns

The frequency

Sets the device list

Parameters

index	the index
id	the ID

```
11.88.3.7 SetImpId() void SetImpId ( uint32\_t id )
```

Sets the ID of the impedance measure device

id the ID


Sets the current RF frequency

Parameters

frequency The frequency

11.89 CPPCDeviceNet Class Reference

Inheritance diagram for CPPCDeviceNet:

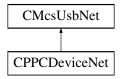

Public Member Functions

• CPPCDeviceNet (void)

Properties

- CPPCFunctionNet^ PPCFunction [get]
- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet^ McsBus_MotorControl [get]
- CMcsBus_SensorNet^ McsBus_Sensor [get]

Additional Inherited Members

11.89.1 Constructor & Destructor Documentation

11.89.1.1 CPPCDeviceNet() CPPCDeviceNet (void)

11.89.2 Property Documentation

```
11.89.2.1 McsBus CMcsBusNet^ McsBus [get]
```

11.89.2.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]


11.89.2.3 McsBus_Sensor CMcsBus_SensorNet^ McsBus_Sensor [get]

11.89.2.4 **PPCFunction** CPPCFunctionNet^ PPCFunction [get]

11.90 CPPCFunctionNet Class Reference

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

Inheritance diagram for CPPCFunctionNet:

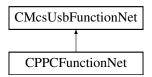

Public Member Functions

CPPCFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pPPCFunctionPointer←
 Container)

Initializes a new instance of the CPPCFunctionNet class.

- CPPCFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CPPCFunctionNet ()
- !CPPCFunctionNet ()
- int GetPumpSpeedUnit (uint16_t channel)

Reads the Pump Speed Unit

void SetPumpSpeedUnit (uint16_t channel, int SpeedUnit)

Writes the Pump Speed Unit

PP_Pump_Mode_Type_EnumNet GetPumpModeType (uint16_t channel)

Reads the Pump Mode Type.

• void SetPumpModeType (uint16_t channel, PP_Pump_Mode_Type_EnumNet PumpMode)

Writes the config string from the device.

void GetAnalogVoltageRange (uint16_t channel, [System::Runtime::InteropServices::Out]uint16_t% min_
 voltage, [System::Runtime::InteropServices::Out]uint16_t% max_voltage)

Reads the Analog Input Voltage Range

void SetAnalogVoltageRange (uint16_t channel, uint16_t min_voltage, uint16_t max_voltage)

Writes the Analog Input Voltage Range

void GetPressureRange (uint16_t channel, [System::Runtime::InteropServices::Out]int32_t% lower_← pressure, [System::Runtime::InteropServices::Out]int32_t% upper pressure)

Get the pressure range that is used between the analog voltage or the digital states

void SetPressureRange (uint16 t channel, int32 t lower pressure, int32 t upper pressure)

Get the pressure range that is used between the analog voltage or the digital states

• uint16_t GetSupplyVoltage ()

Reads the current supply voltage in mV

uint16_t GetAnalogVoltage (uint16_t channel)

Reads the current analog voltage

uint16_t GetDigitalIn (uint16_t channel)

Reads the digital input state

int GetValveActive (uint16_t valve)

Gets the valve active/inactive state

void SetValveActive (uint16 t valve, int valveActive)

Sets the valve active/inactive state

void SetPressureOffset ()

Sets the pressure offset

void LoadPressure (int32 t pressure, uint32 t options)

Loads the reservoir with a pressure

void IsBusy ([System::Runtime::InteropServices::Out]int16_t% task, [System::Runtime::InteropServices::
 Out]int16_t% wait)

Is the PPC busy with a task

• void FirePressurePulse (int32_t duration, int32_t nextpressure)

Fire a pressure pulse from the reservoir

int32_t MeasureReservoir ()

Measures the reservoir pressure

Additional Inherited Members

11.90.1 Detailed Description

CPPCFunctionNet is the class to access the PPC (high precision Patch Peristalic patch Pump

11.90.2 Constructor & Destructor Documentation

Initializes a new instance of the CPPCFunctionNet class.

```
11.90.2.2 CPPCFunctionNet() [2/2] CPPCFunctionNet (
             CMcsUsbNet^ mcsusb )
11.90.2.3 ~CPPCFunctionNet() virtual ~CPPCFunctionNet ( ) [virtual]
11.90.2.4 "!CPPCFunctionNet() !CPPCFunctionNet ( )
11.90.3 Member Function Documentation
11.90.3.1 FirePressurePulse() void FirePressurePulse (
             int32_t duration,
             int32_t nextpressure )
Fire a pressure pulse from the reservoir
Parameters
 duration
               The pulse duration (valves open)
 nextpressure
               The next pressure
11.90.3.2 GetAnalogVoltage() uint16_t GetAnalogVoltage (
             uint16_t channel )
Reads the current analog voltage
Parameters
           The Channel Number
 channel
Returns
     The Analog Voltage
```

11.90.3.3 GetAnalogVoltageRange() void GetAnalogVoltageRange (

uint16_t channel,

```
[System::Runtime::InteropServices::Out] uint16_t% min_voltage, [System::Runtime::InteropServices::Out] uint16_t% max_voltage)
```

Reads the Analog Input Voltage Range

channel	The Channel Number
min_voltage	The voltage that should be seen as the minimum voltage
max_voltage	The voltage that should be seen as the maximum voltage

11.90.3.4 **GetDigitalIn()** uint16_t GetDigitalIn (uint16_t channel)

Reads the digital input state

Parameters

channel The Channel Number

Returns

The Digital State

Get the pressure range that is used between the analog voltage or the digital states

Parameters

channel	The Channel Number
lower_pressure	The lower border of the pressure range
upper_pressure	The upper border of the pressure range

```
11.90.3.6 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType ( uint16_t channel )
```

Reads the Pump Mode Type.

Parameters

channel	The Channel Number

Returns

The Pump Mode Type.

```
11.90.3.7 GetPumpSpeedUnit() int GetPumpSpeedUnit ( uint16_t channel )
```

Reads the Pump Speed Unit

Parameters

channel	The Channel Number
---------	--------------------

Returns

The Speed Unit

11.90.3.8 GetSupplyVoltage() uint16_t GetSupplyVoltage ()

Reads the current supply voltage in mV

Returns

The supply voltage

```
11.90.3.9 GetValveActive() int GetValveActive ( uint16\_t \ valve )
```

Gets the valve active/inactive state

Parameters

```
valve The valve number
```

Returns

The valve state

```
11.90.3.10 IsBusy() void IsBusy (

[System::Runtime::InteropServices::Out] int16_t% task,

[System::Runtime::InteropServices::Out] int16_t% wait )
```

Is the PPC busy with a task

task	The task state
wait	The wait state

Loads the reservoir with a pressure

Parameters

pressure	The pressure
options	The options: end with 0=regulate on patch 1=regulate on reservoir

11.90.3.12 MeasureReservoir() int32_t MeasureReservoir ()

Measures the reservoir pressure

Returns

The pressure

Writes the Analog Input Voltage Range

Parameters

channel	The Channel Number
min_voltage	The voltage that should be seen as the minimum voltage
max_voltage	The voltage that should be seen as the maximum voltage

11.90.3.14 SetPressureOffset() void SetPressureOffset ()

Sets the pressure offset

Get the pressure range that is used between the analog voltage or the digital states

Parameters

channel	The Channel Number
lower_pressure	The lower border of the pressure range
upper_pressure	The upper border of the pressure range

Writes the config string from the device.

Parameters

channel	The Channel Number
PumpMode	The Pump Mode Type.

Writes the Pump Speed Unit

Parameters

channel	The Channel Number
SpeedUnit	The Speed Unit


Sets the valve active/inactive state

Parameters

valve	The valve number
valveActive	The valve state

11.91 CPPS_DeviceNet Class Reference

Inheritance diagram for CPPS_DeviceNet:

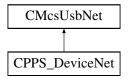

Public Member Functions

CPPS_DeviceNet (void)

Properties

- CPPS_FunctionNet^ PPS_Function [get]
- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet[^] McsBus_MotorControl [get]
- CMcsBus_SensorNet^ McsBus_Sensor [get]

Additional Inherited Members

11.91.1 Constructor & Destructor Documentation

```
11.91.1.1 CPPS_DeviceNet() CPPS_DeviceNet (
void )
```

11.91.2 Property Documentation

```
11.91.2.1 McsBus CMcsBusNet^ McsBus [get]
```


11.91.2.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.91.2.3 McsBus_Sensor CMcsBus_SensorNet^ McsBus_Sensor [get]

11.91.2.4 PPS_Function CPPS_FunctionNet^ PPS_Function [get]

11.92 CPPS FunctionNet Class Reference

Inheritance diagram for CPPS FunctionNet:

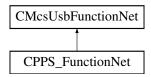

Public Member Functions

- CPPS_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] cPPS_FunctionPointer←
 Container)
- CPPS FunctionNet (CMcsUsbNet[^] mcsusb)
- void SetPumpMaxSpeed (unsigned short index, unsigned short maxspeed)
- unsigned short GetPumpMaxSpeed (unsigned short index)
- void SetPumpSpeedUnit (unsigned short index, int speedunit)
- int GetPumpSpeedUnit (unsigned short index)
- void SetPumpModeType (unsigned short index, PP_Pump_Mode_Type_EnumNet type)
- PP_Pump_Mode_Type_EnumNet GetPumpModeType (unsigned short index)
- void SetPumpCouple (unsigned int i)
- unsigned int GetPumpCouple ()
- void SetPumpEnableSpeedRatio (unsigned int enable)
- unsigned int GetPumpEnableSpeedRatio ()
- void SetPumpManualOnOff (unsigned short index, unsigned int onoff)
- unsigned int GetPumpManualOnOff (unsigned short index)
- · void SetPumpFunctionSpeeds (unsigned short index, short offspeed, short onspeed)
- void GetPumpFunctionSpeeds (unsigned short index, [System::Runtime::InteropServices::Out]short% offspeed, [System::Runtime::InteropServices::Out]short% onspeed)
- void SetPumpSpeedRatio (int ratio)
- int GetPumpSpeedRatio ()
- void SetPumpFastOnOff (unsigned short index, unsigned int onoff)
- unsigned int GetPumpFastOnOff (unsigned short index)
- void SetPumpFastSpeed (unsigned short index, short fastspeed)
- short GetPumpFastSpeed (unsigned short index)
- void SetAnalogVoltages (unsigned short index, unsigned short minvoltage, unsigned short maxvoltage)
- void GetAnalogVoltages (unsigned short index, [System::Runtime::InteropServices::Out]unsigned short% minvoltage, [System::Runtime::InteropServices::Out]unsigned short% maxvoltage)
- void SetUseBubble (unsigned short index, unsigned int usebubble)
- unsigned int GetUseBubble (unsigned short index)
- unsigned short GetSupplyVoltage ()
- unsigned short GetAnalogVoltage (unsigned short index)
- unsigned short GetDigitalIn (unsigned short index)
- unsigned short GetBubbleState ()

Additional Inherited Members

11.92.1 Constructor & Destructor Documentation

```
11.92.1.1 CPPS_FunctionNet() [1/2] CPPS_FunctionNet (
             CMcsUsbNet^ mcsusb,
             {\tt CMcsUsbFunctionPointerContainer}^{\wedge} \ \textit{cPPS\_FunctionPointerContainer} \ )
11.92.1.2 CPPS_FunctionNet() [2/2] CPPS_FunctionNet (
             CMcsUsbNet^ mcsusb )
11.92.2 Member Function Documentation
11.92.2.1 GetAnalogVoltage() unsigned short GetAnalogVoltage (
             unsigned short index )
11.92.2.2 GetAnalogVoltages() void GetAnalogVoltages (
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned short% minvoltage,
             [System::Runtime::InteropServices::Out] unsigned short% maxvoltage )
11.92.2.3 GetBubbleState() unsigned short GetBubbleState ( )
11.92.2.4 GetDigitalIn() unsigned short GetDigitalIn (
             unsigned short index)
11.92.2.5 GetPumpCouple() unsigned int GetPumpCouple ( )
11.92.2.6 GetPumpEnableSpeedRatio() unsigned int GetPumpEnableSpeedRatio ( )
11.92.2.7 GetPumpFastOnOff() unsigned int GetPumpFastOnOff (
             unsigned short index)
```


```
11.92.2.8 GetPumpFastSpeed() short GetPumpFastSpeed (
             unsigned short index)
\textbf{11.92.2.9} \quad \textbf{GetPumpFunctionSpeeds()} \quad \texttt{void GetPumpFunctionSpeeds} \quad \textbf{(}
             unsigned short index,
              [System::Runtime::InteropServices::Out] short% offspeed,
              [System::Runtime::InteropServices::Out] short% onspeed )
11.92.2.10 GetPumpManualOnOff() unsigned int GetPumpManualOnOff (
             unsigned short index)
11.92.2.11 GetPumpMaxSpeed() unsigned short GetPumpMaxSpeed (
             unsigned short index)
11.92.2.12 GetPumpModeType() PP_Pump_Mode_Type_EnumNet GetPumpModeType (
             unsigned short index )
11.92.2.13 GetPumpSpeedRatio() int GetPumpSpeedRatio ( )
11.92.2.14 GetPumpSpeedUnit() int GetPumpSpeedUnit (
             unsigned short index )
11.92.2.15 GetSupplyVoltage() unsigned short GetSupplyVoltage ( )
11.92.2.16 GetUseBubble() unsigned int GetUseBubble (
             unsigned short index )
```

```
11.92.2.17 SetAnalogVoltages() void SetAnalogVoltages (
             unsigned short index,
             unsigned short minvoltage,
             unsigned short {\it maxvoltage} )
11.92.2.18 SetPumpCouple() void SetPumpCouple (
             unsigned int i )
11.92.2.19 SetPumpEnableSpeedRatio() void SetPumpEnableSpeedRatio (
             unsigned int enable )
11.92.2.20 SetPumpFastOnOff() void SetPumpFastOnOff (
             unsigned short index,
             unsigned int onoff )
11.92.2.21 SetPumpFastSpeed() void SetPumpFastSpeed (
             unsigned short index,
             short fastspeed )
11.92.2.22 SetPumpFunctionSpeeds() void SetPumpFunctionSpeeds (
             unsigned short index,
             short offspeed,
             short onspeed )
11.92.2.23 SetPumpManualOnOff() void SetPumpManualOnOff (
             unsigned short index,
             unsigned int onoff )
11.92.2.24 SetPumpMaxSpeed() void SetPumpMaxSpeed (
             unsigned short index,
             unsigned short maxspeed )
```

11.93 CPPSDeviceNet Class Reference

CPPS4plus1DeviceNet is the to control the MCS HLA device

Inheritance diagram for CPPSDeviceNet:

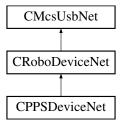

Public Member Functions

• CPPSDeviceNet (void)

Additional Inherited Members

11.93.1 Detailed Description

CPPS4plus1DeviceNet is the to control the MCS HLA device


11.93.2 Constructor & Destructor Documentation

```
11.93.2.1 CPPSDeviceNet() CPPSDeviceNet (
```

11.94 CProgramPressureCurveNet Class Reference

CProgramPressureCurveNet is the class to program pressure curves

Inheritance diagram for CProgramPressureCurveNet:

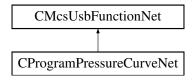

Public Member Functions

- CProgramPressureCurveNet (CMcsUsbNet[^] mcsusb)
 Initializes a new instance of the CPulseGeneratorFunctionNet class.
- virtual ~CProgramPressureCurveNet (void)
- !CProgramPressureCurveNet (void)
- void Program (unsigned char busnumber, unsigned char busaddress, int32_t channel, array< int32_t $>^{\land}$ pressures, array< int32_t $>^{\land}$ durations)
- void SetRepeats (unsigned char busnumber, unsigned char busaddress, int32_t channel, uint32_t repeats)
- unsigned int GetRepeats (unsigned char busnumber, unsigned char busaddress, int32_t channel)

Additional Inherited Members

11.94.1 Detailed Description

CProgramPressureCurveNet is the class to program pressure curves

11.94.2 Constructor & Destructor Documentation

```
11.94.2.1 CProgramPressureCurveNet() CProgramPressureCurveNet ( CMcsUsbNet^ mcsusb )
```


Initializes a new instance of the CPulseGeneratorFunctionNet class.

```
11.94.2.2 ~CProgramPressureCurveNet() virtual ~CProgramPressureCurveNet (
             void ) [virtual]
11.94.2.3 "!CProgramPressureCurveNet() !CProgramPressureCurveNet (
             void )
11.94.3 Member Function Documentation
11.94.3.1 GetRepeats() unsigned int GetRepeats (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel )
11.94.3.2 Program() void Program (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel,
             array< int32_t >^{\wedge} pressures,
             array< int32_t >^{\land} steps,
             array < int16_t >^{\wedge} durations)
11.94.3.3 SetRepeats() void SetRepeats (
             unsigned char busnumber,
             unsigned char busaddress,
             int32_t channel,
             uint32_t repeats )
```

11.95 CPulseGeneratorFunctionNet Class Reference

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

Inheritance diagram for CPulseGeneratorFunctionNet:

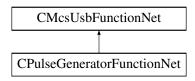

Public Member Functions

 CPulseGeneratorFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pPulse← GeneratorFunctionPointerContainer)

Initializes a new instance of the CPulseGeneratorFunctionNet class.

- CPulseGeneratorFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CPulseGeneratorFunctionNet ()
- !CPulseGeneratorFunctionNet ()
- int32_t GetPeriod (int32_t generator_number)

Reads the generator period

void SetPeriod (int32_t generator_number, int32_t period_in_samples)

Writes the generator period

• int32_t GetPulseLength (int32_t generator_number)

Reads the generator pulse length

void SetPulseLength (int32_t generator_number, int32_t pulselength_in_10us)

Writes the generator pulse length

 void GetModeSelect (int32_t generator_number, [System::Runtime::InteropServices::Out]PulseGenerator_Mode_EnumNet% mode, [System::Runtime::InteropServices::Out]int32_t% digitalchannel)

Reads the generator mode

void SetModeSelect (int32_t generator_number, PulseGenerator_Mode_EnumNet mode, int32_t digitalchannel)

Writes the generator mode

Additional Inherited Members

11.95.1 Detailed Description

CPulseGeneratorFunctionNet is the class to control the pulse generator for video tracking

11.95.2 Constructor & Destructor Documentation

```
11.95.2.1 CPulseGeneratorFunctionNet() [1/2] CPulseGeneratorFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pPulseGeneratorFunctionPointerContainer)
```

Initializes a new instance of the CPulseGeneratorFunctionNet class.

```
11.95.2.2 CPulseGeneratorFunctionNet() [2/2] CPulseGeneratorFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.95.2.3 ~CPulseGeneratorFunctionNet() virtual ~CPulseGeneratorFunctionNet ( ) [virtual]
```

11.95.2.4 "!CPulseGeneratorFunctionNet() !CPulseGeneratorFunctionNet ()

11.95.3 Member Function Documentation

Reads the generator mode

Parameters

generator_number	The generator number
mode	The generator mode
digitalchannel	The digital in channel used as gate

Reads the generator period

Parameters

generator_number	The generator number
------------------	----------------------

Returns

The period

Reads the generator pulse length

Parameters

generator_number	The generator number
------------------	----------------------

Returns

The pulse length

Writes the generator mode

Parameters

generator_number	The generator number
mode	The generator mode
digitalchannel	The digital in channel used as gate

Writes the generator period

Parameters

generator_number	The generator number
period_in_samples	The period


Writes the generator pulse length

Parameters

generator_number	The generator number
pulselength_in_10us	The pulse length

11.96 CRadioControledDevicesNet Class Reference

Inheritance diagram for CRadioControledDevicesNet:

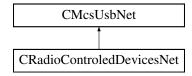

Public Member Functions

- CRadioControledDevicesNet (void)
- bool HasRadioControl ()
- array< unsigned short > ^ GetDeviceNames ()
- void ConnectDevice (unsigned short sn)
- void DisConnectDevice ()
- bool StillConnected ()
- void SetFrequency (unsigned short frequency)
- unsigned short GetFrequency ()

Protected Member Functions

• CRadioControledDevicesNet (CRadioControledDevices *pRadioControled)

Additional Inherited Members

11.96.1 Constructor & Destructor Documentation

```
11.96.1.1 CRadioControledDevicesNet() [1/2] CRadioControledDevicesNet ( void )
```

```
11.96.1.2 CRadioControledDevicesNet() [2/2] CRadioControledDevicesNet (
CRadioControledDevices * pRadioControled) [protected]
```

11.96.2 Member Function Documentation

```
11.96.2.1 ConnectDevice() void ConnectDevice (
unsigned short sn )
```

11.96.2.2 DisConnectDevice() void DisConnectDevice ()

11.97 CCMOSMeaDeviceNet::CRegionOfInterestRect Class Reference

Public Member Functions

- CRegionOfInterestRect (int left, int top, int right, int bottom)
- CRegionOfInterestRect ^ DeepCopy ()

Public Attributes

- int m_Left
- int m_Top
- int m_Right
- int m Bottom

11.97.1 Constructor & Destructor Documentation

```
11.97.1.1 CRegionOfInterestRect() CRegionOfInterestRect (
          int left,
          int top,
          int right,
          int bottom )
```

11.97.2 Member Function Documentation


```
11.97.2.1 DeepCopy() CRegionOfInterestRect ^ DeepCopy ( )
```

11.97.3 Member Data Documentation

```
11.97.3.1 m\_Bottom int m\_Bottom
```

11.98 CRetinaLedDeviceNet Class Reference

Inheritance diagram for CRetinaLedDeviceNet:

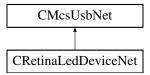

Public Member Functions

- CRetinaLedDeviceNet ()
- ∼CRetinaLedDeviceNet ()
- unsigned int SetTrigger (int enable)
- unsigned int SetLED (unsigned long long pattern)
- unsigned int SetTablepointer (int position)
- unsigned int GetTablepointer (int % position)
- unsigned int ClearTable ()
- unsigned int AddTableEntry (unsigned long long pattern)
- unsigned int AddLoopEntry (unsigned short repeats, unsigned short steps_back)
- unsigned int SetRepeat (int repeat)
- unsigned int SetLumi (int lumi)
- unsigned int SetPersistency (unsigned int persistency)

Additional Inherited Members


```
11.98.1 Constructor & Destructor Documentation
```

```
11.98.1.1 CRetinaLedDeviceNet() CRetinaLedDeviceNet ()
11.98.1.2 ~ CRetinaLedDeviceNet() ~ CRetinaLedDeviceNet ()
11.98.2 Member Function Documentation
11.98.2.1 AddLoopEntry() unsigned int AddLoopEntry (
            unsigned short repeats,
             unsigned short steps_back )
11.98.2.2 AddTableEntry() unsigned int AddTableEntry (
             unsigned long long pattern )
11.98.2.3 ClearTable() unsigned int ClearTable ( )
11.98.2.4 GetTablepointer() unsigned int GetTablepointer (
             int % position )
11.98.2.5 SetLED() unsigned int SetLED (
            unsigned long long pattern )
11.98.2.6 SetLumi() unsigned int SetLumi (
            int lumi )
```

11.99 CRFFunctionNet Class Reference

CRFFunctionNet is the class to control RF devices

Inheritance diagram for CRFFunctionNet:

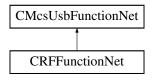

Public Member Functions

CRFFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] pRFFunctionPointer Container

Initializes a new instance of the CRFFunctionNet class.

- CRFFunctionNet (CMcsUsbNet^{\(\Lambda\)} mcsusb)
- virtual ∼CRFFunctionNet ()
- !CRFFunctionNet ()
- uint32_t GetBaseFrequency (CFirmwareDestinationNet destination)

gets the base advertise frequency

void SetBaseFrequency (CFirmwareDestinationNet destination, uint32_t frequency)

sets the base advertise frequency

• uint32_t GetWorkingFrequency ()

gets the working frequency

void SetWorkingFrequency (uint32_t frequency)

sets the working frequency

array< uint32 t > ^ GetAvailableDeviceListEx (int list Length)

get a list of available devices

• array< uint32_t > $^{\wedge}$ GetAvailableDeviceList ()

```
get a list of available devices

    array< uint32_t > ^ GetAvailableStateListEx (int list_Length)

         get a list of the states of the available devices

    array< uint32_t > ^ GetAvailableStateList ()

         get a list of the states of the available devices

    void Connect (uint32_t sn)

         connect to a RF device, use 0 to disconnect
    • uint32_t GetConnectedDevice ()
         get connect RF device, 0 = no device connected

    uint32_t GetState ()

         get connection state

    void SetTestMode (uint32_t mode)

         set test mode
    • uint32_t GetTestMode ()
         gets test mode
Additional Inherited Members
11.99.1 Detailed Description
CRFFunctionNet is the class to control RF devices
11.99.2 Constructor & Destructor Documentation
11.99.2.1 CRFFunctionNet() [1/2] CRFFunctionNet (
               CMcsUsbNet<sup>∧</sup> mcsusb,
               CMcsUsbFunctionPointerContainer^ pRFFunctionPointerContainer )
Initializes a new instance of the CRFFunctionNet class.
11.99.2.2 CRFFunctionNet() [2/2] CRFFunctionNet (
               CMcsUsbNet^ mcsusb )
11.99.2.3 \sim CRFFunctionNet() virtual \sim CRFFunctionNet ( ) [virtual]
11.99.2.4 "!CRFFunctionNet() !CRFFunctionNet ( )
11.99.3 Member Function Documentation
11.99.3.1 Connect() void Connect (
```

uint32_t sn)
connect to a RF device, use 0 to disconnect

Da			_ 1		
Pа	ra	m	eı	re	rs

sn the serial number

```
11.99.3.2 GetAvailableDeviceList() array<uint32_t> ^ GetAvailableDeviceList ( )
```

get a list of available devices

Returns

array of devices

11.99.3.3 GetAvailableDeviceListEx() array<uint32_t> $^{\land}$ GetAvailableDeviceListEx (int $list_Length$)

get a list of available devices

Parameters

list_Length -	The maximal length of list.
---------------	-----------------------------

Returns

array of devices

11.99.3.4 GetAvailableStateList() array<uint32_t> $^{\land}$ GetAvailableStateList ()

get a list of the states of the available devices

Returns

array of states

```
11.99.3.5 GetAvailableStateListEx() array<uint32_t> ^{\land} GetAvailableStateListEx ( int list_Length )
```

get a list of the states of the available devices

list_Length	The maximal length of list.
-------------	-----------------------------

Returns

array of states

```
11.99.3.6 GetBaseFrequency() uint32_t GetBaseFrequency (
CFirmwareDestinationNet destination)
```

gets the base advertise frequency

Parameters

destination the destination to query

Returns

the frequency

11.99.3.7 GetConnectedDevice() uint32_t GetConnectedDevice ()

get connect RF device, 0 = no device connected

Returns

the serial number

11.99.3.8 GetState() uint32_t GetState ()

get connection state

Returns

the state


```
11.99.3.9 GetTestMode() uint32_t GetTestMode ( )
gets test mode
Returns
     the mode
11.99.3.10 GetWorkingFrequency() uint32_t GetWorkingFrequency ( )
gets the working frequency
Returns
     the frequency
11.99.3.11 SetBaseFrequency() void SetBaseFrequency (
             CFirmwareDestinationNet destination,
             uint32_t frequency )
sets the base advertise frequency
Parameters
 destination
             the destination to set
 frequency
             the frequency
11.99.3.12 SetTestMode() void SetTestMode (
             uint32_t mode )
set test mode
Parameters
 mode | the mode
```

```
11.99.3.13 SetWorkingFrequency() void SetWorkingFrequency ( uint32_t frequency )
```

sets the working frequency

11.100 CRobo_FYIProgram_FunctionNet Class Reference

Inheritance diagram for CRobo_FYIProgram_FunctionNet:

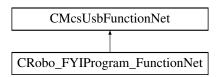

Public Member Functions

- CRobo_FYIProgram_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] robo_←
 FYIProgram FunctionPointerContainer)
- CRobo_FYIProgram_FunctionNet (CMcsUsbNet^ mcsusb)
- void SetValve1 (unsigned char index, unsigned int valve1)
- unsigned int GetValve1 (unsigned char index)
- void SetValve2 (unsigned char index, unsigned int valve2)
- unsigned int GetValve2 (unsigned char index)
- void SetLength (unsigned char index, int length)
- int GetLength (unsigned char index)
- void Start ()
- int GetState ()

Additional Inherited Members

11.100.1 Constructor & Destructor Documentation


```
11.100.1.2 CRobo_FYIProgram_FunctionNet() [2/2] CRobo_FYIProgram_FunctionNet (
CMcsUsbNet^ mcsusb )
```

11.100.2 Member Function Documentation

```
11.100.2.1 GetLength() int GetLength (
             unsigned char index )
11.100.2.2 GetState() int GetState ( )
11.100.2.3 GetValve1() unsigned int GetValve1 (
             unsigned char index )
11.100.2.4 GetValve2() unsigned int GetValve2 (
             unsigned char index)
11.100.2.5 SetLength() void SetLength (
             unsigned char index,
             int length )
11.100.2.6 SetValve1() void SetValve1 (
             unsigned char index,
             unsigned int valve1 )
11.100.2.7 SetValve2() void SetValve2 (
             unsigned char index,
             unsigned int valve2 )
11.100.2.8 Start() void Start ()
```

11.101 CRobo_FYITemp_FunctionNet Class Reference

Inheritance diagram for CRobo_FYITemp_FunctionNet:

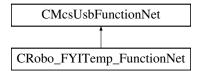

Public Member Functions

- CRobo_FYITemp_FunctionNet (CMcsUsbNet[^] mcsusb)
- void SetRegulatorOnOff (unsigned char index, int onoff)
- int GetRegulatorOnOff (unsigned char index)
- void SetSollTemp (unsigned char index, int temp)
- int GetSollTemp (unsigned char index)
- void SetPCoeff (unsigned char index, int pcoeff)
- int GetPCoeff (unsigned char index)
- void SetlCoeff (unsigned char index, int icoeff)
- int GetlCoeff (unsigned char index)
- void SetMaxPower (unsigned char index, int power)
- int GetMaxPower (unsigned char index)

Additional Inherited Members

11.101.1 Constructor & Destructor Documentation

```
11.101.1.1 CRobo_FYITemp_FunctionNet() CRobo_FYITemp_FunctionNet (
CMcsUsbNet^ mcsusb )
```

11.101.2 Member Function Documentation

```
11.101.2.1 GetlCoeff() int GetlCoeff (
          unsigned char index )
```

```
11.101.2.2 GetMaxPower() int GetMaxPower (
unsigned char index)
```


```
11.101.2.3 GetPCoeff() int GetPCoeff (
unsigned char index )
```

```
11.101.2.4 GetRegulatorOnOff() int GetRegulatorOnOff ( unsigned char index )
```

```
11.101.2.5 GetSollTemp() int GetSollTemp (
             unsigned char index)
11.101.2.6 SetlCoeff() void SetICoeff (
             unsigned char index,
             int icoeff )
11.101.2.7 SetMaxPower() void SetMaxPower (
             unsigned char index,
             int power )
11.101.2.8 SetPCoeff() void SetPCoeff (
             unsigned char index,
             int pcoeff )
11.101.2.9 SetRegulatorOnOff() void SetRegulatorOnOff (
             unsigned char index,
             int onoff )
11.101.2.10 SetSollTemp() void SetSollTemp (
             unsigned char index,
             int temp )
```

11.102 CRoboDacqNet Class Reference

Inheritance diagram for CRoboDacqNet:

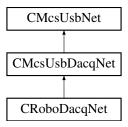

Public Member Functions

- CRoboDacqNet (void)
- CRoboDacqNet (CRoboDeviceNet[^] robodevice)
- void RunTable ()
- void RunTable (int timeout)
- · void StopTable ()
- void StopTable (int timeout)
- void CancelTableLoop ()
- void CancelTableLoopAndStopTable ()
- void SetTriggerMaskValue (unsigned int mask, unsigned int value, unsigned int virtualDevice)
- void SetConfigurationBit (unsigned short bit, bool value)
- void SetConfigurationBitSupply (bool value)
- void SetConfigurationBitRelais (bool value)
- void SetConfigurationBitStream (bool value)
- void SetConfigurationBitAxc (bool value)
- · void SetConfigurationBitCC Gen (bool value)
- void SetConfigurationBitCV_Gen (bool value)
- void SetConfigurationBitRC Gen (bool value)
- void SetConfigurationBitRV_Gen (bool value)
- · void SetConfigurationBitBlu_Led (bool value)
- void SetConfigurationBitRed_Led (bool value)
- void SetConfigurationBitBlu LedToggleSlow (bool value)
- void Octooringulation BitBita_Ecd Toggicolow (bool value)
- void SetConfigurationBitRed_LedToggleSlow (bool value)
- void SetConfigurationBitBlu_LedToggleFast (bool value)
- void SetConfigurationBitRed_LedToggleFast (bool value)void SetConfigurationBitRed_LedSaturation (bool value)
- void SetSimulation (unsigned int enable)
- void SetUClamp (int uClamp)
- void SetIClamp (int iClamp)
- void SetPGain (int pGain)
- void SetIGain (int iGain)
- void SetFilter (int filter)
- void SetUVOffset (int UVOffset)
- void SetUCOffset (int UCOffset)
- void SetICOffset (int ICOffset)
- void SetCrossTalkOffset (int CrossTalk)
- void SetXGain (int xGain)
- void SetCrossTalkOptimum (int cxOptimum)
- void SetRecordingNumber (unsigned int recordingNumber)
- void ClampAmpRestart ()
- void DoRamp (int startValue, int endValue, int duration, int mode)
- unsigned int GetClampAmpSerialNumber ()
- unsigned int GetConfigurationBits ()
- · bool GetConfigurationBit (unsigned short bit)
- bool GetConfigurationBitSupply ()
- bool GetConfigurationBitRelais ()
- bool GetConfigurationBitStream ()
- bool GetConfigurationBitAxc ()
- bool GetConfigurationBitCC_Gen ()
- bool GetConfigurationBitCV_Gen ()
- bool GetConfigurationBitRC Gen ()
- bool GetConfigurationBitRV Gen ()
- bool GetConfigurationBitBlu Led ()
- bool GetConfigurationBitRed_Led ()

- bool GetConfigurationBitBlu_LedToggleSlow ()
- bool GetConfigurationBitRed_LedToggleSlow ()
- bool GetConfigurationBitBlu_LedToggleFast ()
- bool GetConfigurationBitRed_LedToggleFast ()
- bool GetConfigurationBitRed LedSaturation ()
- unsigned int GetSimulation ()
- int GetUClamp ()
- int GetlClamp ()
- int GetPGain ()
- int GetlGain ()
- int GetFilter ()
- int GetUVOffset ()
- int GetUCOffset ()
- int GetICOffset ()
- int dottoonoot ()
- int GetCrossTalkOffset ()
- int GetXGain ()
- int GetCrossTalkOptimum ()
- unsigned int GetRecordingNumber ()
- int GetResistanceC ()
- int GetResistanceV ()
- int GetCapacityC ()
- · int GetCapacityV ()
- · int GetCapacityX ()
- int GetUV ()
- int GetUC ()
- int GetIC ()
- int GetNUV_MS ()
- int GetNUC MS ()
- int GetNIC MS ()
- void SetAllDigout (uint32_t value)
- uint32 t GetAllDigout ()
- void SetCommand (unsigned char command, int value)
- int GetCommand (unsigned char command)
- void SetDigout (uint16 t index, bool enable)
- bool GetDigout (uint16 t index)
- void TableDefBegin ()
- void TableDefEnd ()
- void Table_Wait (unsigned int tableWait)
- void SetDownsampleFactor (int index, int downsample_factor)
- void SetFilterCoeffs (int index, array< int >^ coeffs)
- void SetNoFilterCoeffs (int index)
- int GetDownsampleFactor (int index)
- array< int > ^ GetFilterCoeffs (int index)
- void Emu_SetElectrodeResists (int emuElectrodeResist)
- void Emu_SetCellResists (int emuCellResist)
- void Emu_SetCellCapacity (int emuCellCapacity)
- void Emu SetCellPotential (int emuCellPotential)
- · void Emu_SetNoise (int emuNoise)
- int Emu_GetElectrodeResists ()
- int Emu_GetCellResists ()
- int Emu GetCellCapacity ()
- int Emu_GetCellPotential ()
- int Emu_GetNoise ()
- void SetDisplayText (int index, String[^] displayText)
- void SetScreen (int screen)

- void UpdateDisplay ()
- String \(^\) GetDisplayText (int index)
- int GetScreen ()
- int GetUpdateDisplay ()

Static Public Attributes

- static const unsigned int TriggerMask_Default = 0xFF00
- static const unsigned int TriggerValue_MoveAbs = COMMAND_ROBO_MOVEABS
- static const unsigned int TriggerValue_StartQueue = COMMAND_ROBO_QUEUE
- static const unsigned int VirtualDevice_ContinousDacq = 0
- static const unsigned int VirtualDevice_TableRun = 1

Additional Inherited Members

11.102.1 Constructor & Destructor Documentation

CRoboDeviceNet^ robodevice)

11.102.2 Member Function Documentation

```
11.102.2.1 CancelTableLoop() void CancelTableLoop ( )
```

 $\textbf{11.102.2.2} \quad \textbf{CancelTableLoopAndStopTable()} \quad \texttt{void CancelTableLoopAndStopTable ()} \\$

```
\textbf{11.102.2.3} \quad \textbf{ClampAmpRestart()} \quad \texttt{void ClampAmpRestart ()}
```

```
11.102.2.4 DoRamp() void DoRamp (
             int startValue,
             int endValue,
             int duration,
             int mode )
11.102.2.5 Emu_GetCellCapacity() int Emu_GetCellCapacity ( )
11.102.2.6 Emu_GetCellPotential() int Emu_GetCellPotential ( )
\textbf{11.102.2.7} \quad \textbf{Emu\_GetCellResists()} \quad \texttt{int Emu\_GetCellResists ()}
11.102.2.8 Emu_GetElectrodeResists() int Emu_GetElectrodeResists ( )
11.102.2.9 Emu_GetNoise() int Emu_GetNoise ( )
11.102.2.10 Emu_SetCellCapacity() void Emu_SetCellCapacity (
             int emuCellCapacity )
11.102.2.11 Emu_SetCellPotential() void Emu_SetCellPotential (
             int emuCellPotential )
11.102.2.12 Emu_SetCellResists() void Emu_SetCellResists (
             int emuCellResist )
11.102.2.13 Emu_SetElectrodeResists() void Emu_SetElectrodeResists (
             int emuElectrodeResist )
```

```
11.102.2.14 Emu_SetNoise() void Emu_SetNoise (
            int emuNoise )
11.102.2.15 GetAllDigout() uint32_t GetAllDigout ( )
11.102.2.16 GetCapacityC() int GetCapacityC ( )
11.102.2.17 GetCapacityV() int GetCapacityV ( )
11.102.2.18 GetCapacityX() int GetCapacityX ( )
11.102.2.19 GetClampAmpSerialNumber() unsigned int GetClampAmpSerialNumber ( )
11.102.2.20 GetCommand() int GetCommand (
            unsigned char command )
11.102.2.21 GetConfigurationBit() bool GetConfigurationBit (
            unsigned short bit )
11.102.2.22 GetConfigurationBitAxc() bool GetConfigurationBitAxc ( )
11.102.2.23 GetConfigurationBitBlu_Led() bool GetConfigurationBitBlu_Led ( )
11.102.2.24 GetConfigurationBitBlu_LedToggleFast() bool GetConfigurationBitBlu_LedToggleFast ( )
```

```
11.102.2.25 GetConfigurationBitBlu_LedToggleSlow() bool GetConfigurationBitBlu_LedToggleSlow (
\textbf{11.102.2.26} \quad \textbf{GetConfigurationBitCC\_Gen()} \quad \texttt{bool GetConfigurationBitCC\_Gen ()}
11.102.2.27 GetConfigurationBitCV_Gen() bool GetConfigurationBitCV_Gen ( )
11.102.2.28 GetConfigurationBitRC_Gen() bool GetConfigurationBitRC_Gen ( )
11.102.2.29 GetConfigurationBitRed_Led() bool GetConfigurationBitRed_Led ( )
\textbf{11.102.2.30} \quad \textbf{GetConfigurationBitRed\_LedSaturation()} \quad \texttt{bool GetConfigurationBitRed\_LedSaturation ()} \\
\textbf{11.102.2.31} \quad \textbf{GetConfigurationBitRed\_LedToggleFast()} \quad \texttt{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ( \textbf{bool GetConfigurationBitRed\_LedToggleFast} \ ) \\ \textbf{bool
11.102.2.32 GetConfigurationBitRed_LedToggleSlow() bool GetConfigurationBitRed_LedToggleSlow (
\textbf{11.102.2.33} \quad \textbf{GetConfigurationBitRelais()} \quad \texttt{bool GetConfigurationBitRelais ()}
11.102.2.34 GetConfigurationBitRV_Gen() bool GetConfigurationBitRV_Gen ( )
11.102.2.35 GetConfigurationBits() unsigned int GetConfigurationBits ( )
```

```
11.102.2.36 GetConfigurationBitStream() bool GetConfigurationBitStream ( )
11.102.2.37 GetConfigurationBitSupply() bool GetConfigurationBitSupply ( )
11.102.2.38 GetCrossTalkOffset() int GetCrossTalkOffset ( )
11.102.2.39 GetCrossTalkOptimum() int GetCrossTalkOptimum ( )
11.102.2.40 GetDigout() bool GetDigout (
             uint16_t index )
11.102.2.41 GetDisplayText() String ^ GetDisplayText (
            int index )
11.102.2.42 GetDownsampleFactor() int GetDownsampleFactor (
            int index )
11.102.2.43 GetFilter() int GetFilter ()
11.102.2.44 GetFilterCoeffs() array<int> ^ GetFilterCoeffs (
             int index )
11.102.2.45 GetIC() int GetIC ()
11.102.2.46 GetlClamp() int GetlClamp ( )
```

```
11.102.2.47 GetICOffset() int GetICOffset ( )
11.102.2.48 GetlGain() int GetlGain ()
11.102.2.49 GetNIC_MS() int GetNIC_MS ()
11.102.2.50 GetNUC_MS() int GetNUC_MS ( )
11.102.2.51 GetNUV_MS() int GetNUV_MS ()
11.102.2.52 GetPGain() int GetPGain ( )
\textbf{11.102.2.53} \quad \textbf{GetRecordingNumber()} \quad \texttt{unsigned int GetRecordingNumber ()}
11.102.2.54 GetResistanceC() int GetResistanceC ( )
11.102.2.55 GetResistanceV() int GetResistanceV ()
11.102.2.56 GetScreen() int GetScreen ()
11.102.2.57 GetSimulation() unsigned int GetSimulation ( )
```

```
11.102.2.58 GetUC() int GetUC ()
11.102.2.59 GetUClamp() int GetUClamp ()
11.102.2.60 GetUCOffset() int GetUCOffset ( )
11.102.2.61 GetUpdateDisplay() int GetUpdateDisplay ( )
11.102.2.62 GetUV() int GetUV ( )
11.102.2.63 GetUVOffset() int GetUVOffset ( )
11.102.2.64 GetXGain() int GetXGain ()
11.102.2.65 RunTable() [1/2] void RunTable ( )
11.102.2.66 RunTable() [2/2] void RunTable (
            int timeout )
11.102.2.67 SetAllDigout() void SetAllDigout (
            uint32_t value )
11.102.2.68 SetCommand() void SetCommand (
            unsigned char command,
            int value )
```

```
11.102.2.69 SetConfigurationBit() void SetConfigurationBit (
                                                                                                                                unsigned short bit,
                                                                                                                                bool value )
\textbf{11.102.2.70} \quad \textbf{SetConfigurationBitAxc()} \quad \texttt{void SetConfigurationBitAxc} \quad (
                                                                                                                              bool value )
11.102.2.71 SetConfigurationBitBlu_Led() void SetConfigurationBitBlu_Led (
                                                                                                                                bool value )
\textbf{11.102.2.72} \quad \textbf{SetConfigurationBitBlu\_LedToggleFast()} \quad \texttt{void SetConfigurationBitBlu\_LedToggleFast} \quad \textbf{(} \quad \textbf{(}
                                                                                                                                bool value )
\textbf{11.102.2.73} \quad \textbf{SetConfigurationBitBlu\_LedToggleSlow()} \quad \texttt{void SetConfigurationBitBlu\_LedToggleSlow} \ ( \textbf{void SetConfiguration
                                                                                                                              bool value )
11.102.2.74 SetConfigurationBitCC_Gen() void SetConfigurationBitCC_Gen (
                                                                                                                                bool value )
\textbf{11.102.2.75} \quad \textbf{SetConfigurationBitCV\_Gen()} \quad \texttt{void SetConfigurationBitCV\_Gen()} \quad \texttt{void SetConfigurationBitCV\_Gen()}
                                                                                                                              bool value )
11.102.2.76 SetConfigurationBitRC_Gen() void SetConfigurationBitRC_Gen (
                                                                                                                                bool value )
11.102.2.77 SetConfigurationBitRed_Led() void SetConfigurationBitRed_Led (
                                                                                                                              bool value )
```

```
\textbf{11.102.2.78} \quad \textbf{SetConfigurationBitRed\_LedSaturation()} \quad \texttt{void SetConfigurationBitRed\_LedSaturation} \quad \textbf{(}
                 bool value )
\textbf{11.102.2.79} \quad \textbf{SetConfigurationBitRed\_LedToggleFast()} \quad \texttt{void SetConfigurationBitRed\_LedToggleFast} \ ( \ \textbf{void SetConfigurationBitRed\_LedToggleFast} \ ) \\
                 bool value )
11.102.2.80 SetConfigurationBitRed_LedToggleSlow() void SetConfigurationBitRed_LedToggleSlow (
                bool value )
\textbf{11.102.2.81} \quad \textbf{SetConfigurationBitRelais()} \quad \texttt{void SetConfigurationBitRelais} \quad \textbf{(}
                 bool value )
11.102.2.82 SetConfigurationBitRV_Gen() void SetConfigurationBitRV_Gen (
                 bool value )
11.102.2.83 SetConfigurationBitStream() void SetConfigurationBitStream (
                bool value )
{\bf 11.102.2.84} \quad {\bf SetConfigurationBitSupply()} \quad {\tt void SetConfigurationBitSupply} \ \ (
                 bool value )
11.102.2.85 SetCrossTalkOffset() void SetCrossTalkOffset (
                 int CrossTalk )
11.102.2.86 SetCrossTalkOptimum() void SetCrossTalkOptimum (
                 int cxOptimum )
```

```
11.102.2.87 SetDigout() void SetDigout (
             uint16_t index,
             bool enable )
11.102.2.88 SetDisplayText() void SetDisplayText (
             int index,
             String^{\wedge} displayText)
11.102.2.89 SetDownsampleFactor() void SetDownsampleFactor (
             int index,
             int downsample_factor )
11.102.2.90 SetFilter() void SetFilter (
             int filter )
11.102.2.91 SetFilterCoeffs() void SetFilterCoeffs (
             int index,
             array < int >^{\land} coeffs)
11.102.2.92 SetiClamp() void SetIClamp (
             int iClamp )
11.102.2.93 SetICOffset() void SetICOffset (
             int ICOffset )
11.102.2.94 SetlGain() void SetIGain (
             int iGain )
11.102.2.95 SetNoFilterCoeffs() void SetNoFilterCoeffs (
             int index )
```

```
11.102.2.96 SetPGain() void SetPGain (
            int pGain )
11.102.2.97 SetRecordingNumber() void SetRecordingNumber (
            unsigned int recordingNumber )
11.102.2.98 SetScreen() void SetScreen (
            int screen )
11.102.2.99 SetSimulation() void SetSimulation (
            unsigned int enable )
11.102.2.100 SetTriggerMaskValue() void SetTriggerMaskValue (
             unsigned int mask,
             unsigned int value,
             unsigned int virtualDevice )
11.102.2.101 SetUClamp() void SetUClamp (
            int uClamp )
11.102.2.102 SetUCOffset() void SetUCOffset (
            int UCOffset )
11.102.2.103 SetUVOffset() void SetUVOffset (
            int UVOffset )
11.102.2.104 SetXGain() void SetXGain (
            int xGain )
```


```
11.102.2.105 StopTable() [1/2] void StopTable ( )
11.102.2.106 StopTable() [2/2] void StopTable (
                                    int timeout )
11.102.2.107 Table_Wait() void Table_Wait (
                                    unsigned int tableWait )
11.102.2.108 TableDefBegin() void TableDefBegin ()
11.102.2.109 TableDefEnd() void TableDefEnd ( )
11.102.2.110 UpdateDisplay() void UpdateDisplay ()
11.102.3 Member Data Documentation
11.102.3.1 TriggerMask_Default const unsigned int TriggerMask_Default = 0xFF00 [static]
\textbf{11.102.3.2} \quad \textbf{TriggerValue\_MoveAbs} \quad \texttt{const unsigned int TriggerValue\_MoveAbs} = \texttt{COMMAND\_ROBO\_} \leftarrow
MOVEABS [static]
\textbf{11.102.3.3} \quad \textbf{TriggerValue\_StartQueue} \quad \texttt{const unsigned int TriggerValue\_StartQueue} = \texttt{COMMAND\_ROBO\_} \leftarrow \texttt{COMMAND\_ROBO\_} 
QUEUE [static]
11.102.3.4 VirtualDevice_ContinousDacq const unsigned int VirtualDevice_ContinousDacq = 0
[static]
```

11.102.3.5 VirtualDevice_TableRun const unsigned int VirtualDevice_TableRun = 1 [static]

11.103 CRoboDeviceNet Class Reference

CRoboDeviceNet is the base class for all Robo platform based devices

Inheritance diagram for CRoboDeviceNet:

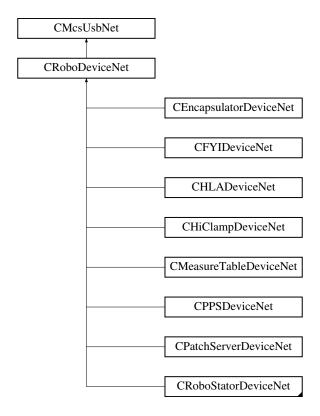

Classes

• class RoboMainLowLevelCommands

Public Member Functions

- CRoboDeviceNet (void)
- ∼CRoboDeviceNet (void)
- void SetInMovement ()

Low level command, sets the internal state to "In Movement"

• bool GetInMovement ()

Low level command, gets the internal state "In Movement"

• uint32_t GetMovementError ()

Low level command, gets the error of the last movement end

- void FindReference (unsigned char busaddress, char axes)
- void FindReference (unsigned char busaddress, char axes, int timeout)

Searches the reference position of the motor

• void MoveAbs (unsigned char busaddress, char axes, int x, int y)

Moves the motor to the new absolute position

void MoveAbs (unsigned char busaddress, char axes, int x, int y, int timeout)

Moves the motor to the new absolute position

void MoveAbs (unsigned char busaddress, char axes, array< int >^ pos)

Moves the motor to the new absolute position

void MoveAbs (unsigned char busaddress, char axes, array< int >^ pos, int timeout)

Moves the motor to the new absolute position

- void StopMovement (unsigned char busaddress, char axes)
- · void StopMovement (unsigned char busaddress, char axes, int timeout)

Stops the current movement

- · void SetCurrentAndAir (unsigned char busaddress, char axes, unsigned short onoff)
- · void SetCurrentAndAir (unsigned char busaddress, char axes, unsigned short onoff, int timeout)
- bool IsQueueEnabled ()
- void EnableQueue (bool enable)
- bool IsQueueStarted ()
- void StartQueue (bool start)
- void WaitTimer (uint32_t waittime, int timeout)
- void CancelPoolLoop ()
- void CancelPoolLoopAndStopMovement ()
- void GetCurrentPosition (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out]int% x, [System::Runtime::InteropServices::Out]int% y)

Gets the current position of motors

- · void SetAirValve (unsigned int onoff)
- unsigned int GetAirValve ()
- void NullCommand (unsigned int marker)
- unsigned int GetVoltageValves ()
- unsigned int GetVoltageRs485A ()
- unsigned int GetVoltageRs485B ()
- unsigned int GetVoltageAirvalve ()
- unsigned int GetCurrentAirvalve ()
- unsigned int GetVoltage12V ()
- unsigned int GetAirpressure ()
- unsigned int GetVoltage5V ()
- unsigned int GetErrorVoltageValves ()
- unsigned int GetErrorVoltageRs485A ()
- unsigned int GetErrorVoltageRs485B ()
- unsigned int GetErrorVoltageAirvalve ()
- unsigned int GetErrorCurrentAirvalve ()
- unsigned int GetErrorVoltage12V ()
- unsigned int GetErrorAirpressure ()
- unsigned int GetErrorVoltage5V ()
- void SetVoltageValvesLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageRs485ALimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageRs485BLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetVoltageAirvalveLimit (unsigned int lowervoltage, unsigned int uppervoltage)
 void SetCurrentAirvalveLimit (unsigned int lowercurrent, unsigned int uppercurrent)
- void SetVoltage12VLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void SetAirpressureLimit (unsigned int lowerpressure, unsigned int upperpressure)
- void SetVoltage5VLimit (unsigned int lowervoltage, unsigned int uppervoltage)
- void GetVoltageValvesLimit ([System::Runtime::InteropServices::Out] unsigned int% lowervoltage, [System ← ::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltageRs485ALimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltageRs485BLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)

- void GetVoltageAirvalveLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetCurrentAirvalveLimit ([System::Runtime::InteropServices::Out]unsigned int% lowercurrent, [System::Runtime::InteropServices::Out]unsigned int% uppercurrent)
- void GetVoltage12VLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::

 Runtime::InteropServices::Out]unsigned int% uppervoltage)
- void GetVoltage5VLimit ([System::Runtime::InteropServices::Out]unsigned int% lowervoltage, [System::
 Runtime::InteropServices::Out]unsigned int% uppervoltage)
- · void SetMinPressure (int pressure)
- int GetMinPressure ()

Static Public Attributes

```
• static const uint32_t RoboError_Base = (0xA0110000L)
```

- static const uint32 t RoboError UnknownCommand = ((0xA0110000L))
- static const uint32 t RoboError Timeout = ((0xA0110000L) | 0x0001)
- static const uint32 t RoboError Pressure = ((0xA0110000L) | 0x0002)
- static const uint32_t RoboError_RangeExceeded = ((0xA0110000L) | 0x0003)
- static const uint32_t RoboError_CommunicationTimeout = ((0xA0110000L) | 0x0004)
- static const uint32_t RoboError_AnotherMaster = ((0xA0110000L) | 0x0005)
- static const uint32 t RoboError FindReferenceMethod = ((0xA0110000L) | 0x0006)
- static const uint32 t RoboError NoSpeedOrAcceleration = ((0xA0110000L) | 0x0007)
- static const uint32_t RoboError_NoEndSwitch = ((0xA0110000L) | 0x0008)
- static const uint32_t RoboError_CannotEscapeEndSwitch = ((0xA0110000L) | 0x0009)
- static const uint32_t RoboError_CommandAlreadyInProgress = ((0xA0110000L) | 0x000A)
- static const uint32_t RoboError_NoReference = ((0xA0110000L) | 0x000B)
- static const uint32_t RoboError_OverPressure = ((0xA0110000L) | 0x000C)
- static const uint32_t RoboError_Phase0OutOfRange = ((0xA0110000L) | 0x000D)
- static const uint32_t RoboError_PeristalticTimeout = ((0xA0110000L) | 0x000E)
- static const uint32_t RoboError_GilsonTimeout = ((0xA0110000L) | 0x000F)
- static const uint32 t RoboError GilsonWrondID = ((0xA0110000L) | 0x0010)
- static const uint32 t RoboError GilsonCommandPending = ((0xA0110000L) | 0x0011)
- static const uint32_t RoboError_ParameterNotAllowed = ((0xA0110000L) | 0x0012)
- static const uint32_t RoboError_StateChangeNotPossible = ((0xA0110000L) | 0x0013)
- static const uint32 t RoboError CommandNotPossible = ((0xA0110000L) | 0x0014)
- static const uint32_t RoboError_DacqNotReady = ((0xA0110000L) | 0x0015)
- static const uint32 t RoboError NoMoreData = ((0xA0110000L) | 0x0016)
- static const uint32_t RoboError_McsBus_UnknownCommand = ((0xA0110000L) | 0x003F)
- static const uint32_t RoboError_DLLMovementTimeout = ((0xA0110000L) | 0x1001)
- static const uint32_t RoboError_PollLoopCanceled = ((0xA0110000L) | 0x1002)
- static const uint32_t RoboError_PollLoopCanceledAndStopMovement = ((0xA0110000L) | 0x1003)
- static const byte McsBus_XY = 1

McsBus address for the xy-plane

• static const byte McsBus ZI = 2

McsBus address for the z and i axes

static const byte Axis_X = 0

Axis number of x for axis argument

static const byte Axis_Y = 1

Axis number of y for axis argument

static const byte Axis_Z = 0

Axis number of z for axis argument

```
• static const byte Axis_I = 1
```

Axis number of i for axis argument

• static const char Axes_X = 1

Bit pattern for x axis for axes argument

static const char Axes_Y = 2

Bit pattern for y axis for axes argument

• static const char Axes Z = 1

Bit pattern for z axis for axes argument

static const char Axes I = 2

Bit pattern for i axis for axes argument

Properties

- CMcsBusNet^ McsBus [get]
- CMcsBus_MotorControlNet^ McsBus_MotorControl [get]
- RoboMainLowLevelCommands A RoboMainLowLevelCommand [get]

Events

• RoboStatusEventDelegate^ RoboStatusEvent

Additional Inherited Members

11.103.1 Detailed Description

CRoboDeviceNet is the base class for all Robo platform based devices

11.103.2 Constructor & Destructor Documentation

```
void )
```

11.103.2.1 CRoboDeviceNet() CRoboDeviceNet (

```
11.103.2.2 ~CRoboDeviceNet() ~CRoboDeviceNet ( void )
```

11.103.3 Member Function Documentation

11.103.3.1 CancelPoolLoop() void CancelPoolLoop ()

```
11.103.3.2 CancelPoolLoopAndStopMovement() void CancelPoolLoopAndStopMovement ( )
```

```
11.103.3.3 EnableQueue() void EnableQueue ( bool enable )
```

```
11.103.3.4 FindReference() [1/2] void FindReference (
unsigned char busaddress,
char axes)
```

Searches the reference position of the motor

Parameters

busaddress	Address of the McsBus	
axes	Bit pattern of axes to drive	
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)	

```
11.103.3.6 GetAirpressure() unsigned int GetAirpressure ( )
```

```
11.103.3.7 GetAirpressureLimit() void GetAirpressureLimit (

[System::Runtime::InteropServices::Out] unsigned int% lowerpressure,

[System::Runtime::InteropServices::Out] unsigned int% upperpressure)
```

11.103.3.8 GetAirValve() unsigned int GetAirValve ()

11.103.3.9 GetCurrentAirvalve() unsigned int GetCurrentAirvalve ()

```
11.103.3.10 GetCurrentAirvalveLimit() void GetCurrentAirvalveLimit (
               [System::Runtime::InteropServices::Out] unsigned int% lowercurrent,
               [{\tt System::Runtime::InteropServices::Out}] \  \, {\tt unsigned int \$} \  \, {\tt uppercurrent} \  \, )
11.103.3.11 GetCurrentPosition() void GetCurrentPosition (
               unsigned char busaddress,
               char axes,
               [System::Runtime::InteropServices::Out] int% x,
               [System::Runtime::InteropServices::Out] int% y )
Gets the current position of motors
Parameters
 busaddress
               Address of the McsBus
                Bit pattern of axes to drive
 axes
                Current position of first axis if pattern in axes is set
 Χ
                Current position of second axis if pattern in axes is set
 у
11.103.3.12 GetErrorAirpressure() unsigned int GetErrorAirpressure ( )
\textbf{11.103.3.13} \quad \textbf{GetErrorCurrentAirvalve()} \quad \texttt{unsigned int GetErrorCurrentAirvalve ()}
\textbf{11.103.3.14} \quad \textbf{GetErrorVoltage12V()} \quad \texttt{unsigned int GetErrorVoltage12V ()}
11.103.3.15 GetErrorVoltage5V() unsigned int GetErrorVoltage5V ()
11.103.3.16 GetErrorVoltageAirvalve() unsigned int GetErrorVoltageAirvalve ()
```

11.103.3.17 GetErrorVoltageRs485A() unsigned int GetErrorVoltageRs485A ()

```
11.103.3.18 GetErrorVoltageRs485B() unsigned int GetErrorVoltageRs485B ( )
11.103.3.19 GetErrorVoltageValves() unsigned int GetErrorVoltageValves ()
11.103.3.20 GetInMovement() bool GetInMovement ( )
Low level command, gets the internal state "In Movement"
11.103.3.21 GetMinPressure() int GetMinPressure ( )
11.103.3.22 GetMovementError() uint32_t GetMovementError ( )
Low level command, gets the error of the last movement end
11.103.3.23 GetVoltage12V() unsigned int GetVoltage12V ()
11.103.3.24 GetVoltage12VLimit() void GetVoltage12VLimit (
              [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
              [System::Runtime::InteropServices::Out] unsigned int% uppervoltage )
11.103.3.25 GetVoltage5V() unsigned int GetVoltage5V ()
11.103.3.26 GetVoltage5VLimit() void GetVoltage5VLimit (
              [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
              [\texttt{System::Runtime::InteropServices::Out}] \ \ unsigned \ \ int \% \ \ uppervoltage \ )
\textbf{11.103.3.27} \quad \textbf{GetVoltageAirvalve()} \quad \texttt{unsigned int GetVoltageAirvalve ()}
```

```
11.103.3.28 GetVoltageAirvalveLimit() void GetVoltageAirvalveLimit (
               [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
               [{\tt System::Runtime::InteropServices::Out}] \ \ unsigned \ \ int {\tt \textit{uppervoltage}} \ )
11.103.3.29 GetVoltageRs485A() unsigned int GetVoltageRs485A ( )
11.103.3.30 GetVoltageRs485ALimit() void GetVoltageRs485ALimit (
               [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
               [\texttt{System::Runtime::InteropServices::Out}] \ \ unsigned \ \ int \% \ \ uppervoltage \ )
11.103.3.31 GetVoltageRs485B() unsigned int GetVoltageRs485B ( )
11.103.3.32 GetVoltageRs485BLimit() void GetVoltageRs485BLimit (
               [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
               [{\tt System::Runtime::InteropServices::Out}] \ \ unsigned \ \ int {\tt \textit{uppervoltage}} \ )
\textbf{11.103.3.33} \quad \textbf{GetVoltageValves()} \quad \texttt{unsigned int GetVoltageValves ()}
11.103.3.34 GetVoltageValvesLimit() void GetVoltageValvesLimit (
               [System::Runtime::InteropServices::Out] unsigned int% lowervoltage,
               [{\tt System::Runtime::InteropServices::Out}] \  \, unsigned \  \, int {\tt \textit{uppervoltage}} \  \, )
11.103.3.35 IsQueueEnabled() bool IsQueueEnabled ( )
11.103.3.36 IsQueueStarted() bool IsQueueStarted ( )
11.103.3.37 MoveAbs() [1/4] void MoveAbs (
               unsigned char busaddress,
               char axes,
              array< int >^{\wedge} pos )
```

Moves the motor to the new absolute position

Parameters

busado	Iress	Address of the McsBus
axes		Bit pattern of axes to drive
pos		Positions of the axis 0 to 3, if pattern in axes is set

```
11.103.3.38 MoveAbs() [2/4] void MoveAbs (
         unsigned char busaddress,
         char axes,
         array< int >^ pos,
         int timeout )
```

Moves the motor to the new absolute position

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
pos	Positions of the axis 0 to 3, if pattern in axes is set
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)

```
11.103.3.39 MoveAbs() [3/4] void MoveAbs (
    unsigned char busaddress,
    char axes,
    int x,
    int y )
```

Moves the motor to the new absolute position

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
X	Position of first axis, if pattern in axes is set
У	Position of second axis if pattern in axes is set

```
11.103.3.40 MoveAbs() [4/4] void MoveAbs (
     unsigned char busaddress,
     char axes,
     int x,
     int y,
     int timeout )
```

Moves the motor to the new absolute position

Parameters

busaddress	Address of the McsBus
axes	Bit pattern of axes to drive
X	Position of first axis, if pattern in axes is set
У	Position of second axis if pattern in axes is set
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)

```
11.103.3.41 NullCommand() void NullCommand (
             unsigned int marker )
11.103.3.42 SetAirpressureLimit() void SetAirpressureLimit (
             unsigned int lowerpressure,
             unsigned int upperpressure )
11.103.3.43 SetAirValve() void SetAirValve (
             unsigned int onoff )
11.103.3.44 SetCurrentAirvalveLimit() void SetCurrentAirvalveLimit (
             unsigned int lowercurrent,
             unsigned int uppercurrent )
11.103.3.45 SetCurrentAndAir() [1/2] void SetCurrentAndAir (
             unsigned char busaddress,
             char axes,
             unsigned short onoff )
11.103.3.46 SetCurrentAndAir() [2/2] void SetCurrentAndAir (
             unsigned char busaddress,
             char axes,
             unsigned short onoff,
             int timeout )
```

```
11.103.3.47 SetInMovement() void SetInMovement ()
Low level command, sets the internal state to "In Movement"
11.103.3.48 SetMinPressure() void SetMinPressure (
              int pressure )
11.103.3.49 SetVoltage12VLimit() void SetVoltage12VLimit (
              unsigned int lowervoltage,
              unsigned int uppervoltage )
11.103.3.50 SetVoltage5VLimit() void SetVoltage5VLimit (
             unsigned int lowervoltage,
              unsigned int uppervoltage )
11.103.3.51 SetVoltageAirvalveLimit() void SetVoltageAirvalveLimit (
              unsigned int lowervoltage,
              unsigned int uppervoltage )
11.103.3.52 SetVoltageRs485ALimit() void SetVoltageRs485ALimit (
              unsigned int lowervoltage,
              unsigned int uppervoltage )
11.103.3.53 SetVoltageRs485BLimit() void SetVoltageRs485BLimit (
              unsigned int lowervoltage,
             unsigned int uppervoltage )
\textbf{11.103.3.54} \quad \textbf{SetVoltageValvesLimit()} \quad \texttt{void SetVoltageValvesLimit ()}
              unsigned int lowervoltage,
```

unsigned int uppervoltage)

```
11.103.3.55 StartQueue() void StartQueue ( bool start )
```

```
11.103.3.56 StopMovement() [1/2] void StopMovement (
unsigned char busaddress,
char axes)
```

```
11.103.3.57 StopMovement() [2/2] void StopMovement (
    unsigned char busaddress,
    char axes,
    int timeout )
```

Stops the current movement

Parameters

busaddress	Address of the McsBus	
axes	Bit pattern of axes to drive	
timeout	Timeout of maximal waiting for the end of the command (-1 is forever)	

11.103.4 Member Data Documentation

```
11.103.4.1 Axes_I const char Axes_I = 2 [static]
```

Bit pattern for i axis for axes argument

```
11.103.4.2 Axes_X const char Axes_X = 1 [static]
```

Bit pattern for x axis for axes argument

```
11.103.4.3 Axes_Y const char Axes_Y = 2 [static]
```

Bit pattern for y axis for axes argument

```
11.103.4.4 Axes_Z const char Axes_Z = 1 [static]
```

Bit pattern for z axis for axes argument

Axis number of i for axis argument

```
11.103.4.6 Axis_X const byte Axis_X = 0 [static]
```

Axis number of x for axis argument

```
11.103.4.7 Axis_Y const byte Axis_Y = 1 [static]
```

Axis number of y for axis argument

```
11.103.4.8 Axis_Z const byte Axis_Z = 0 [static]
```

Axis number of z for axis argument

11.103.4.9 McsBus_XY const byte McsBus_XY = 1 [static]

McsBus address for the xy-plane

11.103.4.10 McsBus_ZI const byte McsBus_ZI = 2 [static]

McsBus address for the z and i axes

```
11.103.4.11 RoboError_AnotherMaster const uint32_t RoboError_AnotherMaster = ( (0xA0110000L) |
0x0005 ) [static]
11.103.4.12 RoboError_Base const uint32_t RoboError_Base = (0xA0110000L) [static]
11.103.4.13 RoboError_CannotEscapeEndSwitch const uint32_t RoboError_CannotEscapeEndSwitch =
( (0xA0110000L) | 0x0009 ) [static]
11.103.4.14 RoboError_CommandAlreadyInProgress const uint32_t RoboError_CommandAlreadyIn←
Progress = ((0xA0110000L) | 0x000A) [static]
11.103.4.15 RoboError_CommandNotPossible const uint32_t RoboError_CommandNotPossible = (
(0xA0110000L) | 0x0014 ) [static]
11.103.4.16 RoboError_CommunicationTimeout const uint32_t RoboError_CommunicationTimeout = (
(0xA0110000L) | 0x0004 ) [static]
11.103.4.17 RoboError DacqNotReady const uint32_t RoboError_DacqNotReady = ( (0xA0110000L) |
0x0015 ) [static]
11.103.4.18 RoboError_DLLMovementTimeout const uint32_t RoboError_DLLMovementTimeout = (
(0xA0110000L) | 0x1001) [static]
11.103.4.19 RoboError_FindReferenceMethod const uint32_t RoboError_FindReferenceMethod = (
(0xA0110000L) | 0x0006 ) [static]
11.103.4.20 RoboError_GilsonCommandPending const uint32_t RoboError_GilsonCommandPending =
( (0xA0110000L) | 0x0011 ) [static]
```

```
11.103.4.21 RoboError_GilsonTimeout const uint32_t RoboError_GilsonTimeout = ( (0xA0110000L) |
0x000F ) [static]
11.103.4.22 RoboError_GilsonWrondID const uint32_t RoboError_GilsonWrondID = ( (0xA0110000L)
| 0x0010 ) [static]
11.103.4.23 RoboError_McsBus_UnknownCommand const uint32_t RoboError_McsBus_Unknown↔
Command = ( (0xA0110000L) | 0x003F) [static]
11.103.4.24 RoboError_NoEndSwitch const uint32_t RoboError_NoEndSwitch = ( (0xA0110000L) |
0x0008 ) [static]
11.103.4.25 RoboError_NoMoreData const uint32_t RoboError_NoMoreData = ( (0xA0110000L) |
0x0016 ) [static]
11.103.4.26 RoboError_NoReference const uint32_t RoboError_NoReference = ( (0xA0110000L) |
0x000B ) [static]
11.103.4.27 RoboError_NoSpeedOrAcceleration const uint32_t RoboError_NoSpeedOrAcceleration =
 ( (0xA0110000L) | 0x0007 ) [static]
11.103.4.28 RoboError_OverPressure const uint32_t RoboError_OverPressure = ( (0xA0110000L) |
0x000C ) [static]
11.103.4.29 RoboError_ParameterNotAllowed const uint32_t RoboError_ParameterNotAllowed = (
 (0xA0110000L) | 0x0012 ) [static]
\textbf{11.103.4.30} \quad \textbf{RoboError\_PeristalticTimeout} \quad \texttt{const uint} \\ 32\_t \quad \texttt{RoboError\_PeristalticTimeout} = \textbf{( } \texttt{(0x} \\ \leftarrow \texttt{)} \\ \textbf{(0x} \\ \leftarrow \texttt{)} \\ \textbf{(0x)} \\ 
A0110000L) | 0x000E ) [static]
```

```
11.103.4.31 RoboError_Phase0OutOfRange const uint32_t RoboError_Phase0OutOfRange = ( (0x↔
A0110000L) | 0x000D ) [static]
11.103.4.32 RoboError_PollLoopCanceled const uint32_t RoboError_PollLoopCanceled = ( (0x \leftarrow
A0110000L) | 0x1002) [static]
\textbf{11.103.4.33} \quad \textbf{RoboError\_PollLoopCanceledAndStopMovement} \quad \texttt{const uint} \\ \textbf{32\_t RoboError\_PollLoop} \leftarrow \\ \textbf{20.1000} \\ \textbf{20.100} \\ 
 CanceledAndStopMovement = ((0xA0110000L) | 0x1003) [static]
11.103.4.34 RoboError_Pressure const uint32_t RoboError_Pressure = ( (0xA0110000L) | 0x0002 )
  [static]
11.103.4.35 RoboError_RangeExceeded const uint32_t RoboError_RangeExceeded = ( (0xA0110000L)
| 0x0003 ) [static]
11.103.4.36 RoboError_StateChangeNotPossible const uint32_t RoboError_StateChangeNotPossible
 = ((0xA0110000L) | 0x0013) [static]
11.103.4.37 RoboError_Timeout const uint32_t RoboError_Timeout = ( (0xA0110000L) | 0x0001 )
  [static]
\textbf{11.103.4.38} \quad \textbf{RoboError\_UnknownCommand} \quad \texttt{const uint} \\ 32\_t \\ \ \textbf{RoboError\_UnknownCommand} \\ = \text{( } \text{(0x} \leftarrow \text{) }
A0110000L) ) [static]
11.103.5 Property Documentation
 11.103.5.1 McsBus CMcsBusNet^ McsBus [get]
```

11.103.5.2 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]


11.103.5.3 RoboMainLowLevelCommand RoboMainLowLevelCommands RoboMainLowLevelCommand [get]

11.103.6 Event Documentation

 $\textbf{11.103.6.1} \quad \textbf{RoboStatusEvent} \quad \texttt{RoboStatusEventDelegate}^{\wedge} \quad \texttt{RoboStatusEvent}$

11.104 CRoboFluidDeviceNet Class Reference

Inheritance diagram for CRoboFluidDeviceNet:

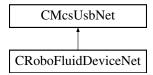

Public Member Functions

- CRoboFluidDeviceNet (void)
- ∼CRoboFluidDeviceNet (void)
- void SetValve (int value)

Open or Close valves.

void SetSingleValve (int valve, bool onoff)

Opens or Closes a valve.

• int GetValve ()

Query the state of the values.

• bool GetSingleValve (int valve)

Query the state of a valve.

- void CloseAllValves ()
- void PumpOn (int index, short speed)
- void SetPumpSpeed (int index, short speed)
- void PumpOff (int index)
- short GetPumpSpeed (int index)
- bool IsPumpMotorOn (int index)

Protected Attributes

- CRoboFluidDevice * m_pRoboFluidDevice
- CMcsBus_MotorControlNet ^ m_pMcsBus_MotorControlNet

Properties

• CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

Additional Inherited Members

11.104.1 Constructor & Destructor Documentation

11.104.2 Member Function Documentation

```
11.104.2.1 CloseAllValves() void CloseAllValves ( )
```

```
11.104.2.2 GetPumpSpeed() short GetPumpSpeed ( int index )
```

```
11.104.2.3 GetSingleValve() bool GetSingleValve ( int valve )
```

Query the state of a valve.

Parameters

valve	number of valve /*!

Returns

state of the valve

```
11.104.2.4 GetValve() int GetValve ()
```

Query the state of the values.

Returns

the current state of the valves as a bit pattern.

```
11.104.2.5 IsPumpMotorOn() bool IsPumpMotorOn ( int index )
```

```
11.104.2.6 PumpOff() void PumpOff ( int index )
```

```
11.104.2.7 PumpOn() void PumpOn (
                int index,
                short speed )
```

Opens or Closes a valve.

Parameters

valve	number of valve to be changed /*!
onoff	open or close the valve

```
11.104.2.10 SetValve() void SetValve ( int value )
```

Open or Close valves.

Parameters

value	bit pattern of valves which should be open.
-------	---

11.104.3 Member Data Documentation

11.104.3.1 m_pMcsBus_MotorControlNet CMcsBus_MotorControlNet ^ m_pMcsBus_MotorControlNet [protected]

11.104.3.2 m_pRoboFluidDevice CRoboFluidDevice* m_pRoboFluidDevice [protected]


11.104.4 Property Documentation

11.104.4.1 McsBus_MotorControl CMcsBus_MotorControlNet^ McsBus_MotorControl [get]

11.105 CRobolnjectDeviceNet Class Reference

CRobolnjectDeviceNet is the to control the MCS Robolnject device

Inheritance diagram for CRobolnjectDeviceNet:

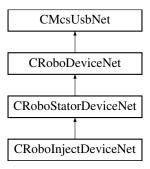

Public Member Functions

CRobolnjectDeviceNet (void)

Additional Inherited Members

11.105.1 Detailed Description

CRobolnjectDeviceNet is the to control the MCS Robolnject device


11.105.2 Constructor & Destructor Documentation

```
11.105.2.1 CRobolnjectDeviceNet() CRobolnjectDeviceNet (
void )
```

11.106 CRoboocyte2DeviceNet Class Reference

CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

Inheritance diagram for CRoboocyte2DeviceNet:

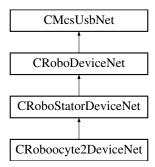

Public Member Functions

- CRoboocyte2DeviceNet (void)
- void SetAxisLED (bool onoff)
- bool GetAxisLED ()
- CRoboDacqNet ^ GetRoboDacq ()
- CRoboFluidDeviceNet ^ GetRoboFluidDevice ()
- CGilsonDeviceNet ^ GetGilsonDevice ()
- CMcsBus_ExtensionNet ^ GetMcsBus_Extension ()

Additional Inherited Members

11.106.1 Detailed Description


CRoboocyte2DeviceNet is the class to control the MCS Roboocyte2 device

11.106.2 Constructor & Destructor Documentation

```
11.106.2.1 CRoboocyte2DeviceNet() CRoboocyte2DeviceNet (
            void )
11.106.3 Member Function Documentation
11.106.3.1 GetAxisLED() bool GetAxisLED ( )
11.106.3.2 GetGilsonDevice() CGilsonDeviceNet ^ GetGilsonDevice ( )
11.106.3.3 GetMcsBus_Extension() CMcsBus_ExtensionNet ^ GetMcsBus_Extension ( )
11.106.3.4 GetRoboDacq() CRoboDacqNet ^ GetRoboDacq ( )
11.106.3.5 GetRoboFluidDevice() CRoboFluidDeviceNet ^ GetRoboFluidDevice ( )
11.106.3.6 SetAxisLED() void SetAxisLED (
            bool onoff )
```

11.107 CRoboStatorDeviceNet Class Reference

Inheritance diagram for CRoboStatorDeviceNet:

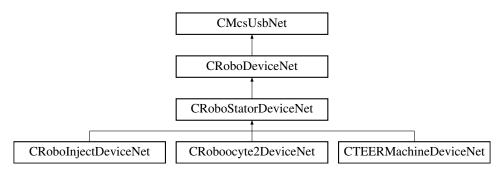

Classes

· class RoboMainStatorLowLevelCommands

Public Member Functions

- CRoboStatorDeviceNet (void)
- void FindReferenceXY ()
- void FindReferenceXY (int timeout)
- void FindReferenceZ ()
- void FindReferenceZ (int timeout)
- void FindReferencel ()
- void FindReferenceI (int timeout)
- unsigned char HasRefXY ()
- unsigned char HasRefZ ()
- unsigned char HasRefl ()
- void MoveAbsXY (int x, int y)
- void MoveAbsXY (int x, int y, int timeout)
- void MoveAbsZ (int z)
- void MoveAbsZ (int z, int timeout)
- void MoveAbsI (int i)
- · void MoveAbsI (int i, int timeout)
- void StopMovementXY ()
- void StopMovementXY (int timeout)
- void StopMovementZ ()
- void StopMovementZ (int timeout)
- void StopMovementI ()
- void StopMovementI (int timeout)
- · void SetCurrentAndAirXY (unsigned short onoff)
- void SetCurrentAndAirXY (unsigned short onoff, int timeout)
- void GetCurrentPositionXY ([System::Runtime::InteropServices::Out]int% x, [System::Runtime::Interop←
 Services::Out]int% y)
- void GetCurrentPositionZ ([System::Runtime::InteropServices::Out]int% z)
- void GetCurrentPositionI ([System::Runtime::InteropServices::Out]int% i)
- void SetVelocityXY (int v)
- void SetVelocityZ (int v)
- void SetVelocityI (int v)
- void SetSpeedXY (int v)
- void SetSpeedZ (int v)
- void SetSpeedI (int v)
- void SetSpeedNativeXY (int v)
- void SetSpeedNativeZ (int v)
- void SetSpeedNativel (int v)
- void SetAccelerationXY (int a)
- void SetAccelerationZ (int a)
- void SetAccelerationI (int a)
- · void SetAccelerationNativeXY (int a)
- void SetAccelerationNativeZ (int a)
- · void SetAccelerationNativel (int a)

Properties

• RoboMainStatorLowLevelCommands RoboMainStatorLowLevelCommand [get]

Additional Inherited Members

11.107.1 Constructor & Destructor Documentation


```
11.107.1.1 CRoboStatorDeviceNet() CRoboStatorDeviceNet (
               void )
11.107.2 Member Function Documentation
11.107.2.1 FindReferencel() [1/2] void FindReferenceI ( )
11.107.2.2 FindReferencel() [2/2] void FindReferenceI (
               int timeout )
\textbf{11.107.2.3} \quad \textbf{FindReferenceXY()} \; \texttt{[1/2]} \quad \texttt{void} \; \texttt{FindReferenceXY} \; \; ( \ )
11.107.2.4 FindReferenceXY() [2/2] void FindReferenceXY (
               int timeout )
11.107.2.5 FindReferenceZ() [1/2] void FindReferenceZ ( )
11.107.2.6 FindReferenceZ() [2/2] void FindReferenceZ (
               int timeout )
\textbf{11.107.2.7} \quad \textbf{GetCurrentPositionI()} \quad \texttt{void GetCurrentPositionI} \quad (
                [System::Runtime::InteropServices::Out] int% i )
```

```
11.107.2.8 GetCurrentPositionXY() void GetCurrentPositionXY (
              [System::Runtime::InteropServices::Out] int% x,
              [System::Runtime::InteropServices::Out] int% y )
\textbf{11.107.2.9} \quad \textbf{GetCurrentPositionZ()} \quad \texttt{void GetCurrentPositionZ} \quad (
              [System::Runtime::InteropServices::Out] int% z )
11.107.2.10 HasRefl() unsigned char HasRefI ()
11.107.2.11 HasRefXY() unsigned char HasRefXY ( )
11.107.2.12 HasRefZ() unsigned char HasRefZ ( )
11.107.2.13 MoveAbsl() [1/2] void MoveAbsI (
             int i)
11.107.2.14 MoveAbsI() [2/2] void MoveAbsI (
             int i,
             int timeout )
11.107.2.15 MoveAbsXY() [1/2] void MoveAbsXY (
             int x,
             int y)
11.107.2.16 MoveAbsXY() [2/2] void MoveAbsXY (
             int x,
             int y,
             int timeout )
```

```
11.107.2.17 MoveAbsZ() [1/2] void MoveAbsZ (
              int z)
11.107.2.18 MoveAbsZ() [2/2] void MoveAbsZ (
             int z,
              int timeout )
11.107.2.19 SetAccelerationI() void SetAccelerationI (
              int a)
11.107.2.20 SetAccelerationNativel() void SetAccelerationNativeI (
              int a)
\textbf{11.107.2.21} \quad \textbf{SetAccelerationNativeXY()} \quad \texttt{void SetAccelerationNativeXY} \quad \textbf{(}
              int a)
11.107.2.22 SetAccelerationNativeZ() void SetAccelerationNativeZ (
              int a )
11.107.2.23 SetAccelerationXY() void SetAccelerationXY (
             int a)
11.107.2.24 SetAccelerationZ() void SetAccelerationZ (
              int a)
11.107.2.25 SetCurrentAndAirXY() [1/2] void SetCurrentAndAirXY (
              unsigned short onoff )
```

```
11.107.2.26 SetCurrentAndAirXY() [2/2] void SetCurrentAndAirXY (
              unsigned short onoff,
              int timeout )
11.107.2.27 SetSpeedI() void SetSpeedI (
             int v )
11.107.2.28 SetSpeedNativel() void SetSpeedNativeI (
             int v)
\textbf{11.107.2.29} \quad \textbf{SetSpeedNativeXY()} \quad \texttt{void SetSpeedNativeXY} \quad \textbf{(}
              int v )
11.107.2.30 SetSpeedNativeZ() void SetSpeedNativeZ (
             int v)
11.107.2.31 SetSpeedXY() void SetSpeedXY (
             int v)
11.107.2.32 SetSpeedZ() void SetSpeedZ (
             int v )
11.107.2.33 SetVelocityI() void SetVelocityI (
              int v)
11.107.2.34 SetVelocityXY() void SetVelocityXY (
             int v)
```

```
11.107.2.35 SetVelocityZ() void SetVelocityZ (
             int v)
11.107.2.36 StopMovementl() [1/2] void StopMovementI ( )
11.107.2.37 StopMovementI() [2/2] void StopMovementI (
             int timeout )
11.107.2.38 StopMovementXY() [1/2] void StopMovementXY ( )
\textbf{11.107.2.39} \quad \textbf{StopMovementXY() [2/2]} \quad \texttt{void StopMovementXY} \quad \textbf{(}
             int timeout )
11.107.2.40 StopMovementZ() [1/2] void StopMovementZ ( )
11.107.2.41 StopMovementZ() [2/2] void StopMovementZ (
             int timeout )
11.107.3 Property Documentation
11.107.3.1 RoboMainStatorLowLevelCommand RoboMainStatorLowLevelCommands ^ RoboMainStator←
LowLevelCommand [get]
11.108 CSafelSDeviceNet Class Reference
Inheritance diagram for CSafeISDeviceNet:
```

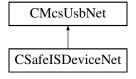

Public Member Functions

CSafeISDeviceNet (void)

Initializes a new instance of the CSafeISDeviceNet class.

∼CSafeISDeviceNet (void)

Releases unmanaged resources and performs other cleanup operations before the CSafelSDeviceNet is reclaimed by garbage collection.

void SetSwitches (unsigned short switches)

Sets the switches for all electrodes on the device. Do not use during measurement

void SetAdcChannels (unsigned char channels)

Sets the ADC channels you want to be sampled

void SetAdcSamplePos (array< unsigned short >[∧] positions)

Sets the sample position of the ADC.

void SetDacMode (unsigned char mode)

Sets the DAC mode.

void SetDacPulseform (array< short >^ pulseform)

Sets the DAC pulseform.

• void SetDacPeriode (unsigned int periode)

Sets the DAC periode.

Properties

• CRoboDeviceNet^ RoboDevice [get]

Gets the CRoboDeviceNet. Use this to control the syringe.

CFluidControlDeviceNet[^] FluidControlDevice [get]

Gets the CFluidControlDeviceNet. Use this to control the valves. Only SetSingleValve is implemented for CSafeISDeviceNet.

• CMcsUsbDacqNet[^] DacqDevice [get]

Gets the CMcsUsbDacqNet. Use this to control the data aquisition.

Additional Inherited Members

11.108.1 Detailed Description

11.108.2 Constructor & Destructor Documentation

```
11.108.2.1 CSafeISDeviceNet() CSafeISDeviceNet (
```

Initializes a new instance of the CSafeISDeviceNet class.

```
11.108.2.2 ~CSafeISDeviceNet() ~CSafeISDeviceNet ( void )
```

Releases unmanaged resources and performs other cleanup operations before the CSafeISDeviceNet is reclaimed by garbage collection.

11.108.3 Member Function Documentation

```
11.108.3.1 SetAdcChannels() void SetAdcChannels (
unsigned char channels )
```

Sets the ADC channels you want to be sampled

Parameters

channels The bitmap of the 8 channels. Set bit to 1 for the channels you want measure

Sets the sample position of the ADC.

Parameters

positions | The positions in units of 0.1μs.

```
11.108.3.3 SetDacMode() void SetDacMode (
unsigned char mode )
```

Sets the DAC mode.

Parameters

mode The mode: 0 = Impedance; 1 = Amperometry

Sets the DAC periode.

Parameters

periode The periode in units of 10μs.

Sets the DAC pulseform.

Parameters

pulseform The pulseform.

```
11.108.3.6 SetSwitches() void SetSwitches (
unsigned short switches)
```

Sets the switches for all electrodes on the device. Do not use during measurement

Parameters

switches The switches: See Schematics for the meaning

11.108.4 Property Documentation

```
11.108.4.1 DacqDevice CMcsUsbDacqNet^ DacqDevice [get]
```

Gets the CMcsUsbDacqNet. Use this to control the data aquisition.

```
11.108.4.2 FluidControlDevice CFluidControlDeviceNet^ FluidControlDevice [get]
```


Gets the CFluidControlDeviceNet. Use this to control the valves. Only SetSingleValve is implemented for CSafeISDeviceNet.

```
11.108.4.3 RoboDevice CRoboDeviceNet^ RoboDevice [get]
```

Gets the CRoboDeviceNet. Use this to control the syringe.

11.109 CSCUDacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CSCUDacqGroupChannelSelectionNet:

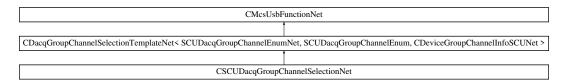

Public Member Functions

CSCUDacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members


11.109.1 Constructor & Destructor Documentation

```
11.109.1.1 CSCUDacqGroupChannelSelectionNet() CSCUDacqGroupChannelSelectionNet (
CMcsUsbNet^ mcsusb)
```

11.110 CSCUFunctionNet Class Reference

CSCUFunctionNet is the class to control the SCU device

Inheritance diagram for CSCUFunctionNet:

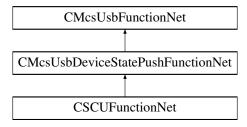

Public Member Functions

- delegate void OnGetAvailableHeadstages (uint32_t AvailableHeadstages)
- delegate void OnlsHeadstageAvailable (uint32_t Headstage, bool available)
- CSCUFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pSCUFunctionPointer←
 Container)

Initializes a new instance of the CSCUFunctionNet class.

- CSCUFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ∼CSCUFunctionNet ()
- !CSCUFunctionNet ()
- uint32 t GetAvailableHeadstages ()

Gets a bitmap of available headstages.

bool IsInDacqLegacyMode ()

Is the SCU in legacy mode

void SetDacqLegacyMode (bool enable)

Enable the SCU legacy mode

uint32_t GetMaxStimulusChannelsPerHeadstage ()

Gets the maximal number of stimulation channels a headstage can have.

uint32 t GetMaxNumberOfHeadstages ()

Gets the maximal number of headstages.

SCU_HeadstageIdEnumNet GetHeadstageID (uint32_t Headstage)

Gets the headstage fpga ID.

bool IsHeadstageAvailable (uint32_t Headstage)

Checks whether the given headstage is available.

void PowerHS (uint32 t Headstage, bool power)

Power the HS

• bool IsHSPowered (uint32 t Headstage)

Is the HS powered

bool HasHSPowerSwitch ()

Has SCU HS power switch

String \(^\) GetHeadstageSerialNumber (uint32_t Headstage)

Gets the serial number of a given headstage.

uint32_t GetHeadstageNumberOfAnalogChannels (uint32_t Headstage)

Gets the number of analog channels for a given headstage.

• void SetHeadstageNumberOfAnalogChannelsPermanent (uint32_t Headstage, uint32_t NumberOfchannels)

Sets the number of analog channels permanent for a given headstage.

uint32_t GetHeadstageNumberOfStimulationChannels (uint32_t Headstage)

Gets the number of stimulation channels for a given headstage.

uint32_t GetHeadstageGainInPermille (uint32_t Headstage)

Gets the gain factor in permille for a given headstage.

uint32 t GetHeadstageAdcRangeInMicroVolt (uint32 t Headstage)

Gets the ADC Range in uV for a given headstage.

uint32_t GetHeadstageAdcBits (uint32_t Headstage)

Gets the Number of ADC bits for a given headstage.

uint32_t GetHeadstageDacVoltageRangeInMilliVolt (uint32_t Headstage)

Gets the DAC Voltage Range in mV for a given headstage.

• uint32_t GetHeadstageDacVoltageResolutionInMicroVolt (uint32_t Headstage)

Gets the DAC Voltage Resolution in uV for a given headstage.

uint32_t GetHeadstageDacCurrentRangeInMicroAmpere (uint32_t Headstage)

Gets the DAC Current Range in uA for a given headstage.

uint32 t GetHeadstageDacCurrentResolutionInNanoAmpere (uint32 t Headstage)

Gets the DAC Current Resolution in nA for a given headstage.

uint32_t GetHeadstageDacBits (uint32_t Headstage)

Gets the Number of DAC bits for a given headstage.

uint32_t GetHeadstageSamplerate (uint32_t Headstage)

Gets the Samplerate of a given headstage.

void SetHeadstageSampleratePermanent (uint32_t Headstage, uint32_t Samplerate)

Sets the samplerate permanent on a given headstage.

uint32_t GetHeadstageLinkSpeed (uint32_t Headstage)

Gets the Link speed of a given headstage.

void SetHeadstageLinkSpeedPermanent (uint32 t Headstage, uint32 t LinkSpeed)

Sets the Link speed permanent on a given headstage.

• uint32 t GetHeadstageFrameCyclesToCompare (uint32 t Headstage)

Gets the frame cycles to compare of a given headstage.

void SetHeadstageFrameCyclesToComparePermanent (uint32_t Headstage, uint32_t FrameCycles)

Sets the frame cycles to compare permanent on a given headstage.

bool GetHeadstagePowerStateAtStart (uint32_t Headstage)

Gets the Power Status at SCU Power on of a given headstage.

void SetHeadstagePowerStateAtStart (uint32_t Headstage, bool Powerstatus)

Sets the Power Status at SCU Power on of a given headstage.

bool HasGalvanicIsolation ()

Has galvanic isolated hardware

bool HasAnalogOut ()

Has AnalogOut hardware

void EnableAnalogOut (bool enable)

Enables AnalogOut globally

bool IsAnalogOutEnabled ()

Is AnalogOut enabled

void SetAnalogOutDACRange (AnalogOut_DAC_Range_EnumNet range)

Sets the analog out DAC range

AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ()

Gets the analog out DAC range

void SetAnalogOutADCRange (uint32_t range)

Sets the analog out ADC range

uint32_t GetAnalogOutADCRange ()

Gets the analog out ADC range

void AutomaticAnalogOut (bool automatic)

Sets automatic source channel selection

bool IsAutomaticAnalogOut ()

Is Automatic source channel selection selected

void SetAnalogOutChannels (uint32_t out_channel, uint32_t source_channel)

Set the source channel number for a certain output channel

uint32_t GetAnalogOutChannels (uint32_t out_channel)

Get the connected source channel number for a certain output channel

 void SetReferenceElectrodeSwitchState (uint32_t Headstage, ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos)

Sets the position of the switch for the reference electrode

• ReferenceElectrodeSwitchPositionEnumNet GetReferenceElectrodeSwitchState (uint32 t Headstage)

Gets the position of the switch for the reference electrode

void SetReferenceElectrodeMode (uint32_t Headstage, ReferenceElectrodeModeEnumNet NewValue)

Sets the mode for the reference electrode

ReferenceElectrodeModeEnumNet GetReferenceElectrodeMode (uint32 t Headstage)

Gets the mode for the reference electrode

- CFilterPropertyNet ^ GetFilterProperty (SCUDacqGroupChannelEnumNet GroupID, uint32_t FilterNumber)
 Gets the filter property
- array< CFilterPropertyNet[^]> [^] GetFilterProperties (SCUDacqGroupChannelEnumNet GroupID, int filter
 — Configurations_Length)

Gets multiple filter properties

Events

• OnGetAvailableHeadstages GetAvailableHeadstagesEvent [add, remove, raise]

Event fires when the bitmap of available headstages has changed

• OnlsHeadstageAvailable^ IsHeadstageAvailableEvent [add, remove, raise]

Event fires when 'true' if the headstage is connected for the headstage to query has changed

Additional Inherited Members

11.110.1 Detailed Description

CSCUFunctionNet is the class to control the SCU device

11.110.2 Constructor & Destructor Documentation

```
11.110.2.1 CSCUFunctionNet() [1/2] CSCUFunctionNet (
             CMcsUsbNet^ mcsusb,
             {\tt CMcsUsbFunctionPointerContainer}^{\wedge} \ p{\tt SCUFunctionPointerContainer} \ )
Initializes a new instance of the CSCUFunctionNet class.
11.110.2.2 CSCUFunctionNet() [2/2] CSCUFunctionNet (
             CMcsUsbNet^ mcsusb )
11.110.2.3 ~CSCUFunctionNet() virtual ~CSCUFunctionNet () [virtual]
11.110.2.4 "!CSCUFunctionNet() !CSCUFunctionNet ()
11.110.3 Member Function Documentation
11.110.3.1 AutomaticAnalogOut() void AutomaticAnalogOut (
             bool automatic )
Sets automatic source channel selection
Parameters
 automatic Automatic
```

Enables AnalogOut globally

Parameters

enable Enable

out_channel Output channel number

Returns

Source channel number

```
11.110.3.5 GetAnalogOutDACRange() AnalogOut_DAC_Range_EnumNet GetAnalogOutDACRange ( )
```

Gets the analog out DAC range

Returns

Range

```
\textbf{11.110.3.6} \quad \textbf{GetAvailableHeadstages()} \quad \texttt{uint32\_t GetAvailableHeadstages ()}
```

Gets a bitmap of available headstages.

Returns

The bitmap of available headstages.

Gets multiple filter properties

GroupID	The group ID
filterConfigurations_Length	The maximal length of filterConfigurations.

Returns

array of filter properties

Gets the filter property

Parameters

GroupID	The group ID
FilterNumber	The filter number

Returns

The filter property

```
11.110.3.9 GetHeadstageAdcBits() uint32_t GetHeadstageAdcBits ( uint32_t Headstage )
```

Gets the Number of ADC bits for a given headstage.

Parameters

Headstage	The headstage to query.

Returns

The number of bits the ADC has for the given headstage.

```
11.110.3.10 GetHeadstageAdcRangeInMicroVolt() uint32_t GetHeadstageAdcRangeInMicroVolt ( uint32_t Headstage )
```

Gets the ADC Range in uV for a given headstage.

Headstage	The headstage to query.
-----------	-------------------------

Returns

The ADC Range in uV for the given headstage.

11.110.3.11 GetHeadstageDacBits() uint32_t GetHeadstageDacBits (uint32_t Headstage)

Gets the Number of DAC bits for a given headstage.

Parameters

Headstage

Returns

The number of bits the DAC has for the given headstage.

11.110.3.12 GetHeadstageDacCurrentRangeInMicroAmpere() uint32_t GetHeadstageDacCurrentRange← InMicroAmpere (uint32_t Headstage)

Gets the DAC Current Range in uA for a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The DAC Current Range in uA for the given headstage.

Gets the DAC Current Resolution in nA for a given headstage.

Returns

The DAC Current Resolution in nA for the given headstage.

Gets the DAC Voltage Range in mV for a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The DAC Voltage Range in mV for the given headstage.

Gets the DAC Voltage Resolution in uV for a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The DAC Voltage Resolution in uV for the given headstage.

11.110.3.16 **GetHeadstageFrameCyclesToCompare()** uint32_t GetHeadstageFrameCyclesToCompare (uint32_t Headstage)

Gets the frame cycles to compare of a given headstage.

Headstage	The headstage to query.
-----------	-------------------------

Returns

The samplerate in Hz for the given headstage.

11.110.3.17 GetHeadstageGainInPermille() uint32_t GetHeadstageGainInPermille (uint32_t Headstage)

Gets the gain factor in permille for a given headstage.

Parameters

Headstage

Returns

The gain factor in permille for the given headstage.

11.110.3.18 **GetHeadstageID()** SCU_HeadstageIdEnumNet GetHeadstageID (uint32_t Headstage)

Gets the headstage fpga ID.

Parameters

Headstage	The headstage to query.

Returns

The headstage fpga ID.

11.110.3.19 **GetHeadstageLinkSpeed()** uint32_t GetHeadstageLinkSpeed (uint32_t Headstage)

Gets the Link speed of a given headstage.

Parameters

Headstage The headstage t

Returns

The samplerate in Hz for the given headstage.

Gets the number of analog channels for a given headstage.

Parameters

Headstage The headstage to que	ry.
--------------------------------	-----

Returns

The number of analog channels the headstage has.

11.110.3.21 GetHeadstageNumberOfStimulationChannels() uint32_t GetHeadstageNumberOfStimulation← Channels (uint32_t Headstage)

Gets the number of stimulation channels for a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The number of stimulation channels the headstage has.

11.110.3.22 GetHeadstagePowerStateAtStart() bool GetHeadstagePowerStateAtStart (uint32_t Headstage)

Gets the Power Status at SCU Power on of a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The Power State at startup for the given headstage: bool false -> off, bool true -> on.

11.110.3.23 GetHeadstageSamplerate() uint32_t GetHeadstageSamplerate (uint32_t Headstage)

Gets the Samplerate of a given headstage.

Parameters

Headstage	The headstage to query.
-----------	-------------------------

Returns

The samplerate in Hz for the given headstage.

11.110.3.24 GetHeadstageSerialNumber() String ^ GetHeadstageSerialNumber (uint32_t Headstage)

Gets the serial number of a given headstage.

Parameters

Headstage The headstage to query.	y.
-----------------------------------	----

Returns

The serial number of the headstage.

11.110.3.25 GetMaxNumberOfHeadstages() uint32_t GetMaxNumberOfHeadstages ()

Gets the maximal number of headstages.

Returns

The maximal number of headstages.

11.110.3.26 GetMaxStimulusChannelsPerHeadstage() uint32_t GetMaxStimulusChannelsPerHeadstage ()

Gets the maximal number of stimulation channels a headstage can have.

Returns

The maximal number of stimulation channels a headstage can have.

Gets the mode for the reference electrode

Parameters

Headstage	The headstage number
-----------	----------------------

Returns

The mode

11.110.3.28 GetReferenceElectrodeSwitchState() ReferenceElectrodeSwitchPositionEnumNet Get↔

Gets the position of the switch for the reference electrode

Parameters

Headstage	The headstage number

Returns

The switch position

11.110.3.29 HasAnalogOut() bool HasAnalogOut ()

Has AnalogOut hardware

Returns

Enabled

```
11.110.3.30 HasGalvanicIsolation() bool HasGalvanicIsolation ()
Has galvanic isolated hardware
Returns
     Enabled
11.110.3.31 HasHSPowerSwitch() bool HasHSPowerSwitch ( )
Has SCU HS power switch
Returns
     Has Switch
11.110.3.32 IsAnalogOutEnabled() bool IsAnalogOutEnabled ( )
Is AnalogOut enabled
Returns
     Enabled
11.110.3.33 IsAutomaticAnalogOut() bool IsAutomaticAnalogOut ( )
Is Automatic source channel selection selected
Returns
     Automatic
11.110.3.34 IsHeadstageAvailable() bool IsHeadstageAvailable (
             uint32_t Headstage )
Checks whether the given headstage is available.
Parameters
              The headstage to query.
 Headstage
```

Returns

'true' if the headstage is connected.

```
11.110.3.35 IsHSPowered() bool IsHSPowered ( uint32_t Headstage )
```

Is the HS powered

Parameters

Returns

'true' if the headstage is powered.

${\tt 11.110.3.36} \quad Is In Dacq Legacy Mode () \quad {\tt bool IsInDacq Legacy Mode ()} \\$

Is the SCU in legacy mode

Returns

Is Enabled

```
11.110.3.37 OnGetAvailableHeadstages() delegate void OnGetAvailableHeadstages ( uint32_t AvailableHeadstages )
```

```
11.110.3.38 OnlsHeadstageAvailable() delegate void OnIsHeadstageAvailable ( uint32_t Headstage, bool available )
```

Power the HS

Headstage	The headstage to query.
power	'true' if the headstage is powered.

11.110.3.40 SetAnalogOutADCRange() void SetAnalogOutADCRange (uint32_t range)

Sets the analog out ADC range

Parameters

range nange		range	Range
-------------	--	-------	-------

Set the source channel number for a certain output channel

Parameters

out_channel	Output channel number
source_channel	Source channel number

Sets the analog out DAC range

Parameters

range	Range

11.110.3.43 SetDacqLegacyMode() void SetDacqLegacyMode (bool *enable*)

Enable the SCU legacy mode

enable	Enable
--------	--------

11.110.3.44 SetHeadstageFrameCyclesToComparePermanent() void SetHeadstageFrameCyclesTo↔

Sets the frame cycles to compare permanent on a given headstage.

Parameters

Headstage	The headstage number
FrameCycles	The samplerate in Hz for the given headstage.

11.110.3.45 SetHeadstageLinkSpeedPermanent() void SetHeadstageLinkSpeedPermanent (

```
uint32_t Headstage,
uint32_t LinkSpeed )
```

Sets the Link speed permanent on a given headstage.

Parameters

Headstage	The headstage number
LinkSpeed	The samplerate in Hz for the given headstage.

$\textbf{11.110.3.46} \quad \textbf{SetHeadstageNumberOfAnalogChannelsPermanent()} \quad \textbf{void SetHeadstageNumberOfAnalog} \leftarrow \textbf{SetHeadstageNumberOfAnalogChannelsPermanent()} \quad \textbf{void SetHeadstageNumberOfAnalogChannelsPermanent()} \quad \textbf{void SetHeadstageNumberOfAnalogChannel$

Sets the number of analog channels permanent for a given headstage.

Parameters

Headstage	The headstage number	
NumberOfchannels	The number of analog channels the headstage has to transmit	l

11.110.3.47 SetHeadstagePowerStateAtStart() void SetHeadstagePowerStateAtStart (

```
uint32_t Headstage,
bool Powerstatus )
```

Sets the Power Status at SCU Power on of a given headstage.

Parameters

Headstage	The headstage number
Powerstatus	The Power State at startup for the given headstage: bool false -> off, bool true -> on.

```
11.110.3.48 SetHeadstageSampleratePermanent() void SetHeadstageSampleratePermanent ( uint32_t Headstage, uint32_t Samplerate)
```

Sets the samplerate permanent on a given headstage.

Parameters

Headstage	The headstage number
Samplerate	The samplerate in Hz for the given headstage.

```
11.110.3.49 SetReferenceElectrodeMode() void SetReferenceElectrodeMode (
    uint32_t Headstage,
    ReferenceElectrodeModeEnumNet NewValue )
```

Sets the mode for the reference electrode

Parameters

Headstage	The headstage number
NewValue	The mode

```
11.110.3.50 SetReferenceElectrodeSwitchState() void SetReferenceElectrodeSwitchState ( uint32_t Headstage,
ReferenceElectrodeSwitchPositionEnumNet NewSwitchPos )
```

Sets the position of the switch for the reference electrode

Parameters

Headstage	The headstage number
NewSwitchPos	The switch position

11.110.4 Event Documentation

11.110.4.1 GetAvailableHeadstagesEvent OnGetAvailableHeadstages^ GetAvailableHeadstagesEvent [add], [remove], [raise]


Event fires when the bitmap of available headstages has changed

11.110.4.2 IsHeadstageAvailableEvent OnIsHeadstageAvailable^ IsHeadstageAvailableEvent [add], [remove], [raise]

Event fires when 'true' if the headstage is connected for the headstage to query has changed

11.111 CSerialPortNet Class Reference

Inheritance diagram for CSerialPortNet:

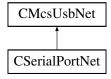

Public Member Functions

- CSerialPortNet (void)
- void Send (array< byte >^ buffer)
- void Send (String[^] command)
- array< byte > $^{\land}$ Receive (void)
- array< byte > ^ Receive (int length)
- String \(^\) ReceiveString (void)
- String \(^\) ReceiveString (int length)
- int GetBytesAvailable (void)

Additional Inherited Members


11.111.1 Constructor & Destructor Documentation

11.111.2 Member Function Documentation

11.112 CStg200xBasicNet Class Reference

Base class for the Stg200x.

Inheritance diagram for CStg200xBasicNet:

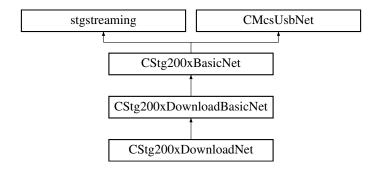

Public Member Functions

virtual ∼CStg200xBasicNet ()

The destructor.

void SetOutputRate (uint32 t rate)

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

uint32_t GetOutputRate ()

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

void SendStart (uint32 t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

• void SendStop (uint32_t triggermap)

Stop some or all triggers of the STG.

void SendStop (uint32 t triggermap, int options)

Stop some or all triggers of the STG.

void GetStgVersionInfo ([Out]String[^]% SwVersion, [Out]String[^]% HwVersion)

Queries software and hardware version.

• void GetAnalogRanges (int channel, [Out]int% URange, [Out]int% IRange)

Gets the range of the analog outputs.

• void GetAnalogResolution (int channel, [Out]int% URes, [Out]int% IRes)

Gets the resolution of the analog outputs.

virtual int32 t GetDACResolution ()

Gets number of bits of the DAC resolution.

virtual int32 t GetTimeResolutionInNanoSeconds ()

Gets the time resolution of the STG.

virtual int32 t GetVoltageRangeInMicroVolt (uint32 t channel)

Gets the Voltage Range of the specified channel in Microvolts.

virtual int32_t GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

virtual int32_t GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

virtual int32_t GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

void GetStgProgramInfo ([Out]bool% IsProgrammed, [Out]System::Runtime::InteropServices::ComTypes::←
FILETIME% timestamp, [Out]String^% filename, [Out]Guid% guid)

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

void GetStgProgramInfo ([Out]bool% IsProgrammed, [Out]DateTime% timestamp, [Out]String^% filename, [Out]Guid% guid)

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

void SetStgProgramInfo (DateTime timestamp, String[^] filename, Guid guid)

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

uint32_t GetAvailableMemory ()

Gets the amount of memory available in the currently selected segment of the STG.

• uint32 t GetTotalMemory ()

Gets the total amount of memory available on the STG (all segments).

virtual uint32_t GetNumberOfAnalogChannels ()

Gets the Number of available analog channels of the device.

virtual uint32 t GetNumberOfSyncoutChannels ()

Gets the Number of available syncout channels of the device.

virtual uint32_t GetNumberOfTriggerInputs ()

Gets the Number of trigger inputs of the device.

virtual uint32_t GetNumberOfHWDACPaths ()

Gets the Number of HW Stimulation DACs of the device.

virtual uint32 t GetNumberOfStimulationSourcesPerElectrode ()

Gets the number of stimulation sources (DACs) per electrode.

virtual void SetVoltageMode (unsigned int channel)

Sets a channel to voltage mode (STG3008-FA and STG400x only).

virtual void SetCurrentMode (unsigned int channel)

Sets a channel to current mode (STG3008-FA and STG400x only).

virtual void SetVoltageMode ()

Sets all channels to voltage mode (STG3008-FA and STG400x only).

• virtual void SetCurrentMode ()

Sets all channels to current mode (STG3008-FA and STG400x only).

virtual System::Collections::Generic::List< int32_t > ^ GetVoltageRangeListInMicroVolt (uint32_t channel)

Gets a list of current ranges supported by the device (STG5 only).

virtual uint32_t GetNumberOfVoltageRangeIndexes ()

Gets the number of voltage ranges (STG5 only).

virtual uint32 t GetVoltageRangeInMicroVoltByIndex (uint32 t channel, uint32 t index)

Gets the voltage range for the given channel and index (STG5 only).

virtual uint32 t GetVoltageResolutionInMicroVoltByIndex (uint32 t channel, uint32 t index)

Gets the voltage resolution for the given channel and index (STG5 only).

virtual uint32_t GetVoltageRangeSelectedIndex (uint32_t channel)

Gets the currently selected range index for the voltage output (not used yet).

virtual void SetVoltageRangeSelectedIndex (uint32_t channel, uint32_t rangeIndex)

Sets the range index for the voltage output (not used yet).

virtual System::Collections::Generic::List< int32 t > ^ GetCurrentRangeListInNanoAmp (uint32 t channel)

Gets a list of current ranges supported by the device (STG5 only).

virtual uint32_t GetNumberOfCurrentRangeIndexes ()

Gets the number of current ranges (STG5 only).

virtual uint32 t GetCurrentRangeInNanoAmpByIndex (uint32 t channel, uint32 t index)

Gets the current range for the given channel and index (STG5 only).

virtual uint32_t GetCurrentResolutionInNanoAmpByIndex (uint32_t channel, uint32_t index)

Gets the current resolution for the given channel and index (STG5 only).

virtual uint32_t GetCurrentRangeSelectedIndex (uint32_t channel)

Gets the currently selected range index for the current output (STG5 only).

virtual void SetCurrentRangeSelectedIndex (uint32_t channel, uint32_t rangeIndex)

Sets the range index for the current output (STG5 only).

virtual void SetMeasurementMode (unsigned int channel)

Sets a channel to measurement mode (STG3008-FA).

virtual void SetFAAmplification (Stg3008FilterAmpAmplificationEnumNet amplification)

Sets the amplification of the filter amplifier (STG3008-FA).

virtual Stg3008FilterAmpAmplificationEnumNet GetFAAmplification ()

Gets the currently selected amplification of the filter amplifier (STG3008-FA).

bool GetCanDisableCalibration (uint32_t channel)

Gets whether the channel can disable autocalibration.

bool GetCanDisableCurrentswitch (uint32_t channel)

Gets whether the channel can disable its current switch.

bool GetHasLowCurrentswitchLatency (uint32_t channel)

Gets whether the channel has a low latency current switch.

virtual void SetAutocalibrationDisabled (unsigned int channel, bool disable)

Sets the autocalibration configuration.

virtual bool GetAutocalibrationDisabled (unsigned int channel)

Gets the autocalibration configuration.

virtual void SetElectrodeMode (uint32_t electrode, array< ElectrodeModeEnumNet >^ mode)

Puts an electrode in either automatic or manual mode.

virtual void SetElectrodeMode (uint32_t electrode, ElectrodeModeEnumNet mode)

Puts an electrode in either automatic or manual mode.

virtual void SetElectrodeMode (uint32_t Scu_HS, uint32_t electrode, array< ElectrodeModeEnumNet >^ mode)

Puts an electrode in either automatic or manual mode.

virtual void SetElectrodeMode (uint32_t Scu_HS, uint32_t electrode, ElectrodeModeEnumNet mode)

Puts an electrode in either automatic or manual mode.

virtual uint32 t GetElectrodeMode (uint32 t electrode)

Gets the mode an electrode is in.

virtual uint32 t GetElectrodeMode (uint32 t Scu HS, uint32 t electrode)

Gets the mode an electrode is in.

virtual void SetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex, array< ElectrodeDacMuxEnumNet
 >^ dacMux)

Defines the DAC to use for an electrode.

virtual void SetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex, ElectrodeDacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

 virtual void SetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, ElectrodeDacMuxEnumNet dacMux)

Defines the DAC to use for an electrode.

virtual void SetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, array
 ElectrodeDacMuxEnumNet >^ dacMux)

Defines the DAC to use for an electrode.

virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t electrode, uint32_t listmodeIndex)

Gets the DAC which is used for an electrode.

virtual ElectrodeDacMuxEnumNet GetElectrodeDacMux (uint32_t Scu_HS, uint32_t electrode, uint32_
 t listmodeIndex)

Gets the DAC which is used for an electrode.

virtual void SetElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an electrode.

virtual void SetElectrodeEnable (uint32_t electrode, uint32_t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an electrode.

virtual void SetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex, array< bool
 <p>>^ enable)

Enables or disables the stimulation switch for an electrode.

• virtual bool GetElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

virtual bool GetElectrodeEnable (uint32_t Scu_HS, uint32_t electrode, uint32_t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

• virtual void SetExternalElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex, array< bool >^ enable)

Enables or disables the stimulation switch for an external electrode.

virtual void SetExternalElectrodeEnable (uint32_t electrode, uint32_t listmodeIndex, bool enable)

Enables or disables the stimulation switch for an external electrode.

virtual bool GetExternalElectrodeEnable (uint32 t electrode, uint32 t listmodeIndex)

Gets weather an electrode is enabled or disabled for stimulation.

virtual void SetBlankingEnable (uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual void SetBlankingEnable (uint32_t electrode, array< bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

• virtual void SetBlankingEnable (uint32_t Scu_HS, uint32_t electrode, bool enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual void SetBlankingEnable (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)

Defines whether an electrode should be blanked while stimulation is in progress.

virtual bool GetBlankingEnable (uint32 t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

virtual bool GetBlankingEnable (uint32_t Scu_HS, uint32_t electrode)

Gets whether an electrode should be blanked while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32_t electrode, bool enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32 t electrode, array < bool >^ enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32 t Scu HS, uint32 t electrode, bool enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual void SetEnableAmplifierProtectionSwitch (uint32_t Scu_HS, uint32_t electrode, array< bool >^ enable)

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual bool GetEnableAmplifierProtectionSwitch (uint32 t electrode)

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

virtual bool GetEnableAmplifierProtectionSwitch (uint32 t Scu HS, uint32 t electrode)

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

- virtual uint32 t GetNumberOfStimulationElectrodes ()
- template<typename digitalsourceenum >

virtual void SetTriggerSource (unsigned int triggernum, DigitalSource< digitalsourceenum >^ triggersource, int bitnum_offset)

- virtual void SetTriggerSource (unsigned int triggernum, TriggerSourceEnumNet triggersource, int bitnum_

 offset)
- virtual void SetTriggerSource (unsigned int triggernum, TriggerSourceEnumNet triggersource)
- virtual TriggerSourceEnumNet GetTriggerSource (unsigned int triggernum)
- virtual void SetListmodeIndexRange (unsigned int electrodeGroup, unsigned int startIndex, unsigned int endIndex, unsigned int mode)

Define the range of list mode indexes to use for the given electrode group.

virtual void GetListmodeIndexRange (unsigned int electrodeGroup, unsigned int &startIndex, unsigned int &endIndex, unsigned int &mode)

Query the range of list mode indexes to use for the given electrode group.

virtual void SetListmodeTriggerSource (unsigned int electrodeGroup, TriggerSourceEnumNet triggersource)

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

virtual void SetListmodeTriggerSource (unsigned int electrodeGroup, TriggerSourceEnumNet triggersource, int bitnumOffset)

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

virtual TriggerSourceEnumNet GetListmodeTriggerSource (unsigned int electrodeGroup)

Query the currently active signal which triggers the transition from one list mode entry to the next. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

virtual void ListModeSendStart (unsigned int electrodeGroupMask)

Activate (arm) the Listmode for the selected electrode groups.

virtual void ListModeSendStop (unsigned int electrodeGroupMask)

Deactivate the Listmode for the selected electrode groups.

- virtual void SetHeadstage (unsigned int headstage)
- virtual uint32_t GetHeadstage ()
- virtual void SetDacAmplificationFactor (uint32 t DacNumber, double Factor)

Set the amplification factor for a DAC.

virtual double GetDacAmplificationFactor (uint32 t DacNumber)

Get the amplification factor for a DAC.

virtual void SetDigoutMode (Stg200xDigoutModeEnumNet digoutMode)

Sets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

virtual Stg200xDigoutModeEnumNet GetDigoutMode()

Gets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

virtual void SetDigoutValue (uint32_t digoutValue)

Sets the Value on the digital output port when in manual mode.

virtual uint32 t GetDigoutValue ()

Gets the Value on the digital output port.

• virtual uint32_t GetDiginValue ()

Gets the Value on the digital input port.

virtual void SetSyncoutMap (uint32_t channel, uint32_t syncoutMap)

Sets the mapping between external syncout outputs and internal syncout channels.

virtual uint32_t GetSyncoutMap (uint32_t channel)

Gets the mapping between external syncout outputs and internal syncout channels.

• virtual bool HasDigitalPort ()

Queries if the device has a digital port.

virtual bool HasFilterAmp ()

Queries if the device has a filter amplifier. Currently the STG3008-FA supports this feature.

virtual bool LegacyHasOutputMap ()

Queries if the device can map its stimulation channels. Legacy command, do not use.

Additional Inherited Members

11.112.1 Detailed Description

Base class for the Stg200x.

From this class all STG related classes are derived: Mcs.Usb.CStg200xDownloadBasicNet Mcs.Usb.CStg200xDownloadNet for Download Mode and Mcs.Usb.CStg200xStreamingNet for Streaming Mode.

CStg200xBasicNet is the base class to control MCS STG device.

11.112.2 Constructor & Destructor Documentation

```
11.112.2.1 \simCStg200xBasicNet() virtual \simCStg200xBasicNet ( ) [virtual]
```

The destructor.

11.112.3 Member Function Documentation

Gets the range of the analog outputs.

Parameters

channel	The channel which is queried.
URange	The Voltage range in mV.
IRange	The Current range in uA.

Gets the resolution of the analog outputs.

Parameters

channel	The channel which is queried.
URes	The Voltage resolution in mV.
IRes	The Current resolution in uA.

```
11.112.3.3 GetAutocalibrationDisabled() virtual bool GetAutocalibrationDisabled ( unsigned int channel ) [virtual]
```

Gets the autocalibration configuration.

Parameters

channel	The channel number.

Returns

true if autocalibration is disabled.

11.112.3.4 GetAvailableMemory() uint32_t GetAvailableMemory ()

Gets the amount of memory available in the currently selected segment of the STG.

Returns

The memory available in the currently selected segment in bytes.

```
11.112.3.5 GetBlankingEnable() [1/2] virtual bool GetBlankingEnable ( uint32_t electrode ) [virtual]
```

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

electrode	The electrode number.
-----------	-----------------------

Returns

true if blanking is enabled while stimulation is in progress.

Gets whether an electrode should be blanked while stimulation is in progress.

Parameters

Scu HS	The SCU headstage number
Scu no	THE SOU HEAUSTAGE HUITIDE

electrode	The electrode number.
electione	I THE EIECHOUE HUITIDEL.

Returns

true if blanking is enabled while stimulation is in progress.

11.112.3.7 GetCanDisableCalibration() bool GetCanDisableCalibration (uint32_t channel)

Gets whether the channel can disable autocalibration.

Parameters

channel	The channel which is queried.
---------	-------------------------------

Returns

True if autocalibration can be disabled.

11.112.3.8 GetCanDisableCurrentswitch() bool GetCanDisableCurrentswitch (uint32_t channel)

Gets whether the channel can disable its current switch.

Parameters

channel	The channel which is queried.

Returns

True if the current switch can be disabled.

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.

Returns

The Current Range of the specified channel in Nanoamps.

Gets the current range for the given channel and index (STG5 only).

Parameters

channel	The channel.
index	The channel.

Returns

The current range in uA.

Gets a list of current ranges supported by the device (STG5 only).

```
11.112.3.12 GetCurrentRangeSelectedIndex() virtual uint32_t GetCurrentRangeSelectedIndex ( uint32_t channel) [virtual]
```

Gets the currently selected range index for the current output (STG5 only).

Parameters

channel	The channel to change.
---------	------------------------

Returns

The currently selected range index.

```
11.112.3.13 GetCurrentResolutionInNanoAmp() virtual int32_t GetCurrentResolutionInNanoAmp ( uint32_t channel ) [virtual]
```

Gets the Current Resolution of the specified channel in Nanoamps.

Returns

The Current Resolution of the specified channel in Nanoamps.

Gets the current resolution for the given channel and index (STG5 only).

Parameters

channel	The channel.
index	The channel.

Returns

The current resolution in nA.

11.112.3.15 GetDacAmplificationFactor() virtual double GetDacAmplificationFactor (uint32_t DacNumber) [virtual]

Get the amplification factor for a DAC.

Parameters

DacNumber	The number of the DAC.
Dadivallibul	

Returns

the amplification factor for the DAC queried, range is from -1.99999 to +1.99999.

11.112.3.16 GetDACResolution() virtual int32_t GetDACResolution () [virtual]

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.112.3.17 GetDiginValue() virtual uint32_t GetDiginValue () [virtual]

Gets the Value on the digital input port.

Returns

The current value on the digital inputs.

11.112.3.18 GetDigoutMode() virtual Stg200xDigoutModeEnumNet GetDigoutMode () [virtual]

Gets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

Returns

The current operation mode.

11.112.3.19 GetDigoutValue() virtual uint32_t GetDigoutValue () [virtual]

Gets the Value on the digital output port.

Returns

The current value on the digital outputs.

Gets the DAC which is used for an electrode.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

Gets the DAC which is used for an electrode.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

The DAC in use, can be 1, 2 or 3. If the electrode is grounded 0 is returned.

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

Gets weather an electrode is enabled or disabled for stimulation.

Scu HS	The SCU headstage number.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

```
11.112.3.24 GetElectrodeMode() [1/2] virtual uint32_t GetElectrodeMode ( uint32_t electrode ) [virtual]
```

Gets the mode an electrode is in.

Parameters

Returns

0 for automatic and 3 for manual mode.

```
11.112.3.25 GetElectrodeMode() [2/2] virtual uint32_t GetElectrodeMode ( uint32_t Scu_HS, uint32_t electrode ) [virtual]
```

Gets the mode an electrode is in.

Scu HS	The SCU headstage number.
000_7.70	i ino oco noddotago nambon

electrode	The electrode number.
-----------	-----------------------

Returns

0 for automatic and 3 for manual mode.

11.112.3.26 GetEnableAmplifierProtectionSwitch() [1/2] virtual bool GetEnableAmplifierProtection← Switch (uint32_t electrode) [virtual]

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

electrode	The electrode number.
-----------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

Gets whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Scu_HS	The SCU headstage number.

electrode	The electrode number.
-----------	-----------------------

Returns

true if the switch is to be opened, false if it is closed while stimulation is in progress.

Gets weather an electrode is enabled or disabled for stimulation.

Parameters

electrode	The electrode number.
listmodeIndex	The index for listmode.

Returns

true if the electrode is enabled, false if it is disabled.

```
11.112.3.29 GetFAAmplification() virtual Stg3008FilterAmpAmplificationEnumNet GetFAAmplification ( ) [virtual]
```

Gets the currently selected amplification of the filter amplifier (STG3008-FA).

Returns

True the currently selected amplification.

```
11.112.3.30 GetHasLowCurrentswitchLatency() bool GetHasLowCurrentswitchLatency ( uint32_t channel)
```

Gets whether the channel has a low latency current switch.

Parameters

channel The channel which is queried.

Returns

True if the channel has a low latency current switch.

```
11.112.3.31 GetHeadstage() virtual uint32_t GetHeadstage ( ) [virtual]
```

11.112.3.32 **GetListmodeIndexRange()** virtual void GetListmodeIndexRange (unsigned int *electrodeGroup*,

```
unsigned int & startIndex,
unsigned int & endIndex,
unsigned int & mode ) [virtual]
```

Query the range of list mode indexes to use for the given electrode group.

Parameters

electrodeGroup	The electrodegroup for which the range is queried.
----------------	--

Parameters

startIndex	The index of the first active element in the listmode list.
endIndex	The index of the last active element in the listmode list.
mode	0 for "start with startIndex", 1 for "start with endIndex".

11.112.3.33 GetListmodeTriggerSource() virtual TriggerSourceEnumNet GetListmodeTriggerSource (unsigned int electrodeGroup) [virtual]

Query the currently active signal which triggers the transition from one list mode entry to the next. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

Parameters

electrodeGroup	The electrodegroup for which the triggersource is queried.
----------------	--

Returns

One of the possible sources for the transition.

11.112.3.34 GetNumberOfAnalogChannels() virtual uint32_t GetNumberOfAnalogChannels () [virtual]

Gets the Number of available analog channels of the device.

Returns

The number of analog channels.

11.112.3.35 GetNumberOfCurrentRangeIndexes() virtual uint32_t GetNumberOfCurrentRangeIndexes () [virtual]

Gets the number of current ranges (STG5 only).

Returns

The number of current ranges available on the device.

11.112.3.36 GetNumberOfHWDACPaths() virtual uint32_t GetNumberOfHWDACPaths () [virtual]

Gets the Number of HW Stimulation DACs of the device.

Returns

The number of independent HW Stimulation outputs.

11.112.3.37 GetNumberOfStimulationElectrodes() virtual uint32_t GetNumberOfStimulationElectrodes () [virtual]

11.112.3.38 GetNumberOfStimulationSourcesPerElectrode() virtual uint32_t GetNumberOfStimulation← SourcesPerElectrode () [virtual]

Gets the number of stimulation sources (DACs) per electrode.

Returns

The number of stimulation sources (DACs) per electrode.

11.112.3.39 **GetNumberOfSyncoutChannels()** virtual uint32_t GetNumberOfSyncoutChannels () [virtual]

Gets the Number of available syncout channels of the device.

Returns

The number of analog channels.

11.112.3.40 GetNumberOfTriggerInputs() virtual uint32_t GetNumberOfTriggerInputs () [virtual]

Gets the Number of trigger inputs of the device.

Returns

The number of trigger inputs.

11.112.3.41 GetNumberOfVoltageRangeIndexes() virtual uint32_t GetNumberOfVoltageRangeIndexes () [virtual]

Gets the number of voltage ranges (STG5 only).

Returns

The number of voltage ranges available on the device.

11.112.3.42 GetOutputRate() uint32_t GetOutputRate ()

Queries the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Returns

Returns the current output rate in Hz.

11.112.3.43 GetStgProgramInfo() [1/2] void GetStgProgramInfo (

```
[Out] bool% IsProgrammed,
[Out] DateTime% timestamp,
[Out] String^% filename,
[Out] Guid% guid)
```

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

IsProgrammed	Flag wether download information is valid.
timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.
guid	A GUID.

11.112.3.44 GetStgProgramInfo() [2/2] void GetStgProgramInfo (

```
[Out] bool% IsProgrammed,
[Out] System::Runtime::InteropServices::ComTypes::FILETIME% timestamp,
[Out] String^% filename,
[Out] Guid% guid )
```

Queries Download information from the STG. If download information was stored by the use of SetStgProgramInfo, this function can be used to retrieve it.

Parameters

IsProgrammed	Flag wether download information is valid.
timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.

Queries software and hardware version.

Parameters

SwVersion	The current Software Version of the STG.
HwVersion	The Hardware Revision of the STG.

Gets the mapping between external syncout outputs and internal syncout channels.

Parameters

channel	The external syncout output channel number (zero based).
---------	--

Returns

The bitmap of internal syncout channels mapped to channel.

```
11.112.3.47 GetTimeResolutionInNanoSeconds() virtual int32_t GetTimeResolutionInNanoSeconds (
) [virtual]
```

Gets the time resolution of the STG.

Returns

The time resolution in nanoseconds.

11.112.3.48 GetTotalMemory() uint32_t GetTotalMemory ()

Gets the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

```
11.112.3.49 GetTriggerSource() virtual TriggerSourceEnumNet GetTriggerSource ( unsigned int triggernum ) [virtual]
```

```
11.112.3.50 GetVoltageRangeInMicroVolt() virtual int32_t GetVoltageRangeInMicroVolt ( uint32_t channel) [virtual]
```

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

channel	Channel which is queried.

Returns

The Voltage Range of the specified channel in Microvolts.

Gets the voltage range for the given channel and index (STG5 only).

channel	The channel.
index	The channel.

Returns

The voltage range.

Gets a list of current ranges supported by the device (STG5 only).

```
11.112.3.53 GetVoltageRangeSelectedIndex() virtual uint32_t GetVoltageRangeSelectedIndex ( uint32_t channel) [virtual]
```

Gets the currently selected range index for the voltage output (not used yet).

Parameters

channel	The channel to change.
---------	------------------------

Returns

The currently selected range index.

```
11.112.3.54 GetVoltageResolutionInMicroVolt() virtual int32_t GetVoltageResolutionInMicroVolt ( uint32_t channel ) [virtual]
```

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

_		
	channel	Channel which is queried.

Returns

The Voltage Resolution of the specified channel in Microvolts.

Gets the voltage resolution for the given channel and index (STG5 only).

channel	The channel.
index	The channel.

Returns

The voltage resolution.

11.112.3.56 HasDigitalPort() virtual bool HasDigitalPort () [virtual]

Queries if the device has a digital port.

Returns

true if the device has a digital port.

11.112.3.57 HasFilterAmp() virtual bool HasFilterAmp () [virtual]

Queries if the device has a filter amplifier. Currently the STG3008-FA supports this feature.

Returns

true if the device has a a filter amplifier.

11.112.3.58 LegacyHasOutputMap() virtual bool LegacyHasOutputMap () [virtual]

Queries if the device can map its stimulation channels. Legacy command, do not use.

Returns

true if the device can map its stimulation channels.

```
 \begin{array}{llll} \textbf{11.112.3.59} & \textbf{ListModeSendStart()} & \texttt{virtual void ListModeSendStart (} \\ & & \texttt{unsigned int } electrodeGroupMask \ ) & \texttt{[virtual]} \end{array}
```

Activate (arm) the Listmode for the selected electrode groups.

electrodeGroupMask	The bitmask of electrode groups for which the listmode is activated.

```
11.112.3.60 ListModeSendStop() virtual void ListModeSendStop (
unsigned int electrodeGroupMask) [virtual]
```

Deactivate the Listmode for the selected electrode groups.

Parameters

ectrodeGroupMask The bitmask of electrodegroup	ps for which the listmode is deactivated.
--	---

```
11.112.3.61 SendStart() void SendStart ( uint32_t triggermap )
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

triggermap	A bitmap of triggers which will be started.
------------	---

```
11.112.3.62 SendStop() [1/2] void SendStop ( uint32_t triggermap )
```

Stop some or all triggers of the STG.

Parameters

triggermap A bitmap of triggers which will be stopped.

```
11.112.3.63 SendStop() [2/2] void SendStop (
    uint32_t triggermap,
    int options)
```

Stop some or all triggers of the STG.

triggermap	A bitmap of triggers which will be stopped.	
options	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses	
	the stop commands when a syncout assossiated with a given sync-out has bit 1 (0x02) set. Can	
	be used e.g. to prevent a stop while a biphasic stimulation pulse is active	

```
11.112.3.64 SetAutocalibrationDisabled() virtual void SetAutocalibrationDisabled (
unsigned int channel,
bool disable) [virtual]
```

Sets the autocalibration configuration.

Parameters

channel	The channel number.
disable	true if autocalibration is to be disabled.

```
11.112.3.65 SetBlankingEnable() [1/4] virtual void SetBlankingEnable ( uint32_t electrode, array< bool >^ enable ) [virtual]
```

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

```
11.112.3.66 SetBlankingEnable() [2/4] virtual void SetBlankingEnable ( uint32_t electrode, bool enable ) [virtual]
```

Defines whether an electrode should be blanked while stimulation is in progress.

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

Scu_HS The SCU headstag	ge number.
-------------------------	------------

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether an electrode should be blanked while stimulation is in progress.

Parameters

Scu_HS	The SCU headstage number.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

11.112.3.69 SetCurrentMode() [1/2] virtual void SetCurrentMode () [virtual]

Sets all channels to current mode (STG3008-FA and STG400x only).

```
11.112.3.70 SetCurrentMode() [2/2] virtual void SetCurrentMode (
unsigned int channel) [virtual]
```

Sets a channel to current mode (STG3008-FA and STG400x only).

The channel to change.	channel
------------------------	---------

11.112.3.71 SetCurrentRangeSelectedIndex() virtual void SetCurrentRangeSelectedIndex (uint32_t channel, uint32_t rangeIndex) [virtual]

Sets the range index for the current output (STG5 only).

Parameters

channel	The channel to change.
rangeIndex	The new range index.

Set the amplification factor for a DAC.

Parameters

DacNumber	The number of the DAC.
Factor	the amplification factor for that DAC, range is from -1.99999 to +1.99999.

```
11.112.3.73 SetDigoutMode() virtual void SetDigoutMode (
Stg200xDigoutModeEnumNet digoutMode) [virtual]
```

Sets the operation mode of the digital outport port, can be Monitor, Manual or SyncOut

Parameters

digoutMode	The new operation mode.
------------	-------------------------

```
11.112.3.74 SetDigoutValue() virtual void SetDigoutValue ( uint32_t digoutValue ) [virtual]
```

Sets the Value on the digital output port when in manual mode.

digoutValue	The new value on the digital outputs.
-------------	---------------------------------------

Defines the DAC to use for an electrode.

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use
	ElectrodeDacMuxEnumNet.Ground (0).

Defines the DAC to use for an electrode.

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

```
\textbf{11.112.3.77} \quad \textbf{SetElectrodeDacMux() [3/4]} \quad \text{virtual void SetElectrodeDacMux (}
```

```
uint32_t Scu_HS,
uint32_t electrode,
uint32_t listmodeIndex,
array< ElectrodeDacMuxEnumNet >^ dacMux ) [virtual]
```

Defines the DAC to use for an electrode.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2 (2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use ElectrodeDacMuxEnumNet.Ground (0).

```
\textbf{11.112.3.78} \quad \textbf{SetElectrodeDacMux() [4/4]} \quad \text{virtual void SetElectrodeDacMux (}
```

```
uint32_t Scu_HS,
uint32_t electrode,
uint32_t listmodeIndex,
ElectrodeDacMuxEnumNet dacMux ) [virtual]
```

Defines the DAC to use for an electrode.

Scu_HS	The SCU headstage number.
--------	---------------------------

electrode The electrode number.

Parameters

listmodeIndex	The index for listmode.
dacMux	The DAC to use, can be ElectrodeDacMuxEnumNet.Stg1 (1), ElectrodeDacMuxEnumNet.Stg2
	(2) or ElectrodeDacMuxEnumNet.Stg3 (3). To ground an electrode, use
	ElectrodeDacMuxEnumNet.Ground (0).

Enables or disables the stimulation switch for an electrode.

Parameters

electrode The electrode nu

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an electrode.

electrode T	he electrode number.
-------------	----------------------

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an electrode.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
- electrode	THE Electrode Hulliber.

Parameters

listmodeInde	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

Enables or disables the stimulation switch for an electrode.

Scu HS	The SCU headstage number.

Parameters

electrode	The electrode number.
-----------	-----------------------

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

```
11.112.3.83 SetElectrodeMode() [1/4] virtual void SetElectrodeMode ( uint32_t electrode, array< ElectrodeModeEnumNet >^ mode ) [virtual]
```

Puts an electrode in either automatic or manual mode.

Parameters

electrode	The electrode number.
electione i	THE ELECTIONS HUTTINGS.

Returns

0 for automatic and 3 for manual mode.

Puts an electrode in either automatic or manual mode.

electrode	The electrode number.
electrode	i The electrode number.

mode 0 for automatic and 3 for man	ual mode.
------------------------------------	-----------

Puts an electrode in either automatic or manual mode.

Parameters

Scu_HS The SCU headstage number.

Parameters

Returns

0 for automatic and 3 for manual mode.

Puts an electrode in either automatic or manual mode.

Scu HS	The SCU headstage number.

electrode	The electrode number.	
-----------	-----------------------	--

Parameters

mode 0 for a	automatic and 3 for manual mode.
--------------	----------------------------------

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

electrode	The electrode number.	
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.	

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

electrode	The electrode number.
enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Defines whether the Amplifier Protection Switch is openend while stimulation is in progress.

Parameters

Scu_HS	The SCU headstage number.
--------	---------------------------

Parameters

Ī	electrode	The electrode number.
	enable	True if the switch is to be opened, false if it is to remain closed while stimulation is in progress.

Enables or disables the stimulation switch for an external electrode.

electrode The electrode number	er.
--------------------------------	-----

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

```
11.112.3.92 SetExternalElectrodeEnable() [2/2] virtual void SetExternalElectrodeEnable (
    uint32_t electrode,
    uint32_t listmodeIndex,
    bool enable ) [virtual]
```

Enables or disables the stimulation switch for an external electrode.

Parameters

electrode The electrode numbe	r.	
-------------------------------	----	--

Parameters

listmodeIndex	The index for listmode.
enable	1 to enable the electrode, 0 to disable.

```
11.112.3.93 SetFAAmplification() virtual void SetFAAmplification (
Stg3008FilterAmpAmplificationEnumNet amplification ) [virtual]
```

Sets the amplification of the filter amplifier (STG3008-FA).

```
11.112.3.94 SetHeadstage() virtual void SetHeadstage (
unsigned int headstage ) [virtual]
```

Define the range of list mode indexes to use for the given electrode group.

<i>electrodeGroup</i> The electrodegroup for which the range is defined.
--

Parameters

startIndex	The index of the first active element in the listmode list.
endIndex	The index of the last active element in the listmode list.
mode	0 for "start with startIndex", 1 for "start with endIndex".

```
11.112.3.96 SetListmodeTriggerSource() [1/2] virtual void SetListmodeTriggerSource (
unsigned int electrodeGroup,
TriggerSourceEnumNet triggersource) [virtual]
```

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

Parameters

electrodeGroup	The electrodegroup for which the triggersource is defined.
cicolicacaioap	The diediedegroup for which the thiggerodated is defined.

Parameters

ti	riggersource	One of the possible sources for the transition.
----	--------------	---

Define the signal which triggers the transition from one list mode entry to the next. After reaching the last entry in the list, the first entry is selected. For triggersource use the Enum which corresponds to the device in use, for example use SCUDigitalSourceEnumNet and cast to TriggerSourceEnumNet if working with an SCU device.

electrodeGroup	The electrodegroup for which the triggersource is defined.

triggersource	One of the possible sources for the transition.	
bitnumOffset	Number to add to the numeric value of the <i>triggersource</i> enum.	

11.112.3.98 SetMeasurementMode() virtual void SetMeasurementMode (unsigned int channel) [virtual]

Sets a channel to measurement mode (STG3008-FA).

Parameters

channel	The channel to change.
---------	------------------------

```
11.112.3.99 SetOutputRate() void SetOutputRate ( uint32_t rate )
```

Change the output rate of the STG. Valid rates are from 1000 Hz to 50000 Hz.

Parameters

rate	The new output rate in Hz.

```
11.112.3.100 SetStgProgramInfo() void SetStgProgramInfo (
DateTime timestamp,
String^ filename,
Guid guid )
```

Store Download information in the STG. This function can be used to store the filename and timestamp of the last download for later query.

Parameters

timestamp	The timestamp of last download.
filename	The filename of the downlaoded waveform.

11.112.3.101 SetSyncoutMap() virtual void SetSyncoutMap (

```
uint32_t channel,
uint32_t syncoutMap ) [virtual]
```

Sets the mapping between external syncout outputs and internal syncout channels.

Parameters

channel	The external syncout output channel number (zero based).
syncoutMap	A bitmap of internal syncout channels to map to channel.

```
11.112.3.102 SetTriggerSource() [1/3] virtual void SetTriggerSource (
          unsigned int triggernum,
          DigitalSource< digitalsourceenum >^ triggersource,
          int bitnum_offset ) [virtual]
```

```
11.112.3.103 SetTriggerSource() [2/3] virtual void SetTriggerSource (
unsigned int triggernum,
TriggerSourceEnumNet triggersource) [virtual]
```

```
11.112.3.104 SetTriggerSource() [3/3] virtual void SetTriggerSource (
    unsigned int triggernum,
    TriggerSourceEnumNet triggersource,
    int bitnum_offset ) [virtual]
```

```
11.112.3.105 SetVoltageMode() [1/2] virtual void SetVoltageMode ( ) [virtual]
```

Sets all channels to voltage mode (STG3008-FA and STG400x only).

```
11.112.3.106 SetVoltageMode() [2/2] virtual void SetVoltageMode ( unsigned int channel ) [virtual]
```

Sets a channel to voltage mode (STG3008-FA and STG400x only).

channel	The channel to change.
---------	------------------------

Sets the range index for the voltage output (not used yet).


Parameters

channel	The channel to change.
rangeIndex	The new range index.

11.113 CStg200xDownloadBasicNet Class Reference

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

Inheritance diagram for CStg200xDownloadBasicNet:

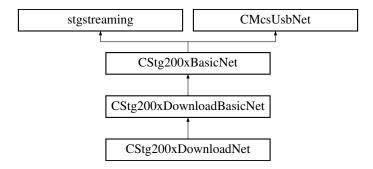

Public Member Functions

virtual void SetupTrigger (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

- virtual void SetupTriggerSingle (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)

 Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.
- void GetTrigger ([Out] array< uint32_t >^% channelmap, [Out] array< uint32_t >^% syncoutmap, [Out] array< uint32_t >^% repeat)

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

- void GetSweepCount ([Out] array< uint32_t >^% sweeps, [Out] array< uint32_t >^% triggers)

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in CStg200xDownloadBasicNet::SetupTrigger.
- void ForceStatusEvent ()

Force a status event.

void ResetStatus (uint32_t triggermap)

Reset the status flag.

uint32 t GetMemoryUsageDAC (uint32 t Channel)

Queries the memory usage of the current segment and selected analog DAC channel.

uint32_t GetMemoryUsageSyncout (uint32_t Channel)

Queries the memory usage of the current segment and selected syncout channel.

virtual void ClearSyncData (uint32_t channel)

Delete a SyncOut pattern for a channel from STG memory.

virtual void SendSyncData (uint32_t channel, array< uint16_t >^ pData, array< uint64_t >^ tData)

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the STG200x_SendSyncData32() call instead.

• virtual void ClearChannelData (uint32 t channel)

Delete a stimulus pattern for a channel from STG memory

• virtual void SendChannelData (uint32_t channel, array< uint16_t $>^{\land}$ pData, array< uint64_t $>^{\land}$ tData)

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μs. The STG has a resolution of 20 μs.

virtual void EnableAutoReset ()

Enable AutoReset of the STG Status.

virtual void DisableAutoReset ()

Disable AutoReset of the STG Status.

virtual void SetupRetriggerMode (int8_t trigger, RetriggerActionEnumNet same_trigger, RetriggerActionEnumNet other trigger)

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal
- virtual void SetupRetriggerMode (RetriggerActionEnumNet same_trigger, RetriggerActionEnumNet other_
 trigger)

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- restart this trigger
- ignore the signal

Properties

CStimulusFunctionNet[^] Stimulus [get]

Additional Inherited Members

11.113.1 Detailed Description

CStg200xDownloadBasicNet is the base class to control the download mode of the MCS STG device.

11.113.2 Member Function Documentation

11.113.2.1 ClearChannelData() virtual void ClearChannelData (uint32_t channel) [virtual]

Delete a stimulus pattern for a channel from STG memory

Parameters

channel Specifies the channel to clear.

11.113.2.2 ClearSyncData() virtual void ClearSyncData (uint32_t channel) [virtual]

Delete a SyncOut pattern for a channel from STG memory.

Parameters

channel Specifies the syncout channel to clear.

11.113.2.3 DisableAutoReset() virtual void DisableAutoReset () [virtual]

Disable AutoReset of the STG Status.

If autoreset is disabled, the STG status switches to FINISHED after the defined number of sweeps is finished. To switch back to the IDLE status, use CStg200xDownload::ResetStatus()

11.113.2.4 EnableAutoReset() virtual void EnableAutoReset () [virtual]

Enable AutoReset of the STG Status.

This is the default on power up. If autoreset is enabled, the STG status switches to FINISHED only for one poll cycle after this, it switches to IDLE automatically.

11.113.2.5 ForceStatusEvent() void ForceStatusEvent ()

Force a status event.

Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

Queries the memory usage of the current segment and selected analog DAC channel.

The currently used memory is reported for the requested channel.

Channel	channel for the amount of interested usage.
---------	---

Returns

Returns the usage in STG memory.

```
11.113.2.7 GetMemoryUsageSyncout() uint32_t GetMemoryUsageSyncout ( uint32_t Channel)
```

Queries the memory usage of the current segment and selected syncout channel.

The currently used memory is reported for the requested channel.

Parameters

	Channel	channel for the amount of interested usage.	
--	---------	---	--

Returns

Returns the usage in STG memory.

```
11.113.2.8 GetSweepCount() void GetSweepCount (

[Out] array< uint32_t >^% sweeps,

[Out] array< uint32_t >^% triggers )
```

Get the sweep and trigger count of the STG.

- The triggercount tells how many times each trigger was active and is reset to zero on download of new channel data.
- The sweepcount tells how many times each trigger was already repeated. This count is set to zero on trigger and counts up to repeat in CStg200xDownloadBasicNet::SetupTrigger.

sweeps	on return contains the number of sweeps for each trigger.
triggers	on return contains the number of trigger events seen for each trigger.

```
11.113.2.9 GetTrigger() void GetTrigger (

[Out] array< uint32_t >^% channelmap,
```

```
[Out] array< uint32_t >^{\%} syncoutmap, [Out] array< uint32_t >^{\%} repeat )
```

Queries the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

channelmap	For each trigger, a bitmap of channels that belong to this trigger.
------------	---

Parameters

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

11.113.2.10 ResetStatus() void ResetStatus (uint32_t triggermap)

Reset the status flag.

Parameters

```
triggermap bitmap of trigger for which to reset the status.
```

Uploads analog data (stimulus patterns) to the STG.

Sends datapoints to a given channel on the STG. The list of datapoints will be sent to the selected channel. Data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value in the range from 0 to 4095 (bit 0 to 11), its sign is taken from bit 12, 0 is for positive amplitude, and 1 for negative amplitude Bits 13 to 15 have to be zero.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

channel Specifies the channel to append the data to.	
pData	A list of datapoints.
tData	A list of durations as int64_t. The time is given in units of μs.

Uploads sync output data to the STG.

Sends sync output data to a given channel on the STG. The list of datapoints will be sent to the selected sync output channel. Sync output data previously sent to the channel is overwritten.

Each datapoint is represented by an integer value and can be either 0 or 1.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s. If your application can not handle 64 bit integers, use the STG200x_SendSyncData32() call instead.

Parameters

channel	Specifies the sync output channel to append the data to.
pData	A list of datapoints.
tData	A list of durations as int64_t. The time is given in units of μs.

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- · restart this trigger
- · ignore the signal

Parameters

trigger	The trigger to change.
same_trigger	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
other_trigger	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected. Not used in Normal Mode.

11.113.2.14 SetupRetriggerMode() [2/2] virtual void SetupRetriggerMode (

```
RetriggerActionEnumNet same_trigger,
RetriggerActionEnumNet other_trigger) [virtual]
```

Define the action on triggers while the STG is running.

The STG has three options how to handle a successive trigger while a trigger is active.

- stop this trigger (default action)
- · restart this trigger
- · ignore the signal

Parameters

same_trigger	Action for successive triggers in Normal Mode, and for triggers to the currently selected segment in Multi - File Mode.
other_trigger	Action for successive triggers in Multi-File Mode for a trigger on a segment not currently selected.Not used in Normal Mode.

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

	first_trigger	The number of the first trigger to change.
--	---------------	--

Parameters

channelmap	For each trigger, a bitmap of channels that belong to this trigger.
------------	---

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

11.113.2.16 SetupTriggerSingle() virtual void SetupTriggerSingle (

```
uint32_t trigger,
uint32_t channelmap,
uint32_t syncoutmap,
uint32_t repeat ) [virtual]
```

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

trigger	The trigger to change.
---------	------------------------

Parameters

Parameters

syncoutmap	A bitmap of syncouts that belong to this trigger.
repeat	The number of times this trigger should be repeated.


11.113.3 Property Documentation

```
11.113.3.1 Stimulus CStimulusFunctionNet^ Stimulus [get]
```

11.114 CStg200xDownloadNet Class Reference

Main class for the STG download mode This class implements the STG download mode interface.

Inheritance diagram for CStg200xDownloadNet:

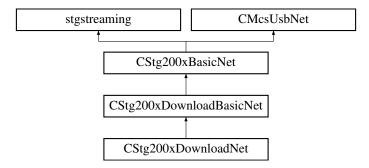

Public Member Functions

CStg200xDownloadNet ()

Use this constructor if you do not want to use the status callback.

CStg200xDownloadNet (OnStgPollStatus[^] pollStatus)

Use this constructor if you want to use the status callback.

- ∼CStg200xDownloadNet ()
- void PrepareAndSendData (uint32_t channel, array< int32_t $>^{\land}$ amplitude, array< uint64_t $>^{\land}$ duration, STG DestinationEnumNet destType)

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

 void PrepareAndAppendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and append data to a given channel on the STG.

void ClearChannel_PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType, bool doClear)

Prepare and append data to a given channel on the STG.

void SegmentDefine (array< uint32_t >^ segment_list)

Defines the segment memory layout of the STG.

- void SegmentStart (uint32_t triggermap, uint32_t segment, Stg200xSegmentFlagsEnumNet segmentflags)
 Switchs segment and starts trigger.
- void SegmentSelect (uint32 t segment, Stg200xSegmentFlagsEnumNet segmentflags)

Switchs segment.

• void EnableMultiFileMode (Stg200xMultiFileSubmodeEnumNet submode)

Enable the Multi-File mode of the STG.

· void DisableMultiFileMode ()

Disable the Multi-File mode of the STG

- StgStatusNet ^ QueryTriggerstatus ()
- void SetOutputMap (array< uint32_t >^ ChannelLayout)
- int32_t GetModuleTemp (unsigned int channel)
- uint32_t GetModuleCurrent (unsigned int channel)

Events

- OnStgPollStatus^ Stg200xPollStatusEvent [add, remove, raise]
- OnMwPollStatus^ MwPollStatusEvent [add, remove, raise]

Additional Inherited Members

11.114.1 Detailed Description

Main class for the STG download mode This class implements the STG download mode interface.

11.114.2 Constructor & Destructor Documentation

11.114.2.1 CStg200xDownloadNet() [1/2] CStg200xDownloadNet ()

Use this constructor if you do not want to use the status callback.

```
11.114.2.2 CStg200xDownloadNet() [2/2] CStg200xDownloadNet (
OnStgPollStatus^ pollStatus)
```

Use this constructor if you want to use the status callback.

```
11.114.2.3 ~CStg200xDownloadNet() ~CStg200xDownloadNet ()
```

11.114.3 Member Function Documentation

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

end data to.	number to	The channel	channel
--------------	-----------	-------------	---------

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

duration	A list of durations in units of μs .
destType	specifies wheather the data is for syncout, current or voltage stimulation.

11.114.3.2 DisableMultiFileMode() void DisableMultiFileMode ()

Disable the Multi-File mode of the STG

Switch the STG back to normal mode. In this mode, trigger inputs are assigned to channels, not to segments.

```
11.114.3.3 EnableMultiFileMode() void EnableMultiFileMode ( Stg200xMultiFileSubmodeEnumNet submode )
```

Enable the Multi-File mode of the STG.

In Multi-File mode, the trigger inputs switch between segments. To use this mode, define up to as many segments as trigger inputs are available and fill each segment with a stimulus pattern.

Now a trigger on trigger input 1 switches the STG to the first segment and starts all triggers in this segment. Likewise, a trigger on trigger input 2, 3 and 4 selects the respective segment and start all triggers in this segment So the Multi-File Mode can be used to predefine up to four different stimuli which can be selected without the need for a computer connection.

submode	The submode. Submode 0 is regular Multi-File mode as described above, submode 1 is extended
	Multi-File mode, where the segment is selected based on the digital pattern on the digital inputs. In
	this mode, 256 different segments can be defined and used.

```
11.114.3.5 GetModuleTemp() int32_t GetModuleTemp ( unsigned int channel )
```

11.114.3.6 PrepareAndAppendData() void PrepareAndAppendData (

```
uint32_t channel,
array< int32_t >^ amplitude,
array< uint64_t >^ duration,
STG_DestinationEnumNet destType )
```

Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel The channel number to se	end data to.
----------------------------------	--------------

Parameters

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

Parameters

duration	A list of durations in units of μs.	1
destType	specifies wheather the data is for syncout, current or voltage stimulation.	

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel	The channel number to send data to.
---------	-------------------------------------

Parameters

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

Parameters

duration	A list of durations in units of μs.
destType	specifies wheather the data is for syncout, current or voltage stimulation.

11.114.3.8 QueryTriggerstatus() StgStatusNet ^ QueryTriggerstatus ()

```
11.114.3.9 SegmentDefine() void SegmentDefine ( array < uint32_t >^{\land} segment_list )
```

Defines the segment memory layout of the STG.

On reset, the STG has one segment containing all available memory.

With this command, the STG memory can be devided into several segments. Each segment can be filled with stimulus data.

```
segment_list  The List of memory sizes (one per segment).
```

Switchs segment.

Parameters

	segment	The number of the segment to select.
--	---------	--------------------------------------

Parameters

segmentflags A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.

Switchs segment and starts trigger.

Parameters

triggermap	A bitmap of triggers that will be started.
------------	--

Parameters

segment	The number of the segment to select.
---------	--------------------------------------

Parameters

segmentflags A bitmap of flags, bit 1: assign all channels to the trigger number equal to the segment.

```
11.114.3.12 SetOutputMap() void SetOutputMap (

array< uint32_t >^ ChannelLayout )
```


11.114.4 Event Documentation

11.114.4.1 MwPollStatusEvent OnMwPollStatus^ MwPollStatusEvent [add], [remove], [raise]

11.114.4.2 Stg200xPollStatusEvent OnStgPollStatus^ Stg200xPollStatusEvent [add], [remove], [raise]

11.115 CStimulusFunctionNet Class Reference

Inheritance diagram for CStimulusFunctionNet:

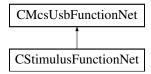

Classes

- · class SidebandData
- · class StimulusDeviceDataAndUnrolledData

Public Member Functions

- CStimulusFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] stimulusFunction PointerContainer)
- CStimulusFunctionNet (CMcsUsbNet[∧] mcsusb)
- void StartPoll ()

Starts the interrupt fetching thread and delivers events

void StopPoll ()

Stops the interrupt fetching thread and delivers events

void ForceStatusEvent ()

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

void SendStart (uint32 t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

void SendStop (uint32_t triggermap)

Stop some or all triggers of the STG.

void SendStop (uint32_t triggermap, int options)

Stop some or all triggers of the STG.

void ClearChannelData (int channel)

Delete a Stimulus Pattern from STG memory

void ClearSyncData (int channel)

Delete a Syncout Pattern from STG memory

• void PrepareAndSendData (uint32_t channel, array< int32_t $>^{\land}$ amplitude, array< uint64_t $>^{\land}$ duration, STG_DestinationEnumNet destType)

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

 void PrepareAndAppendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType)

Prepare and append data to a given channel on the STG.

- void ClearChannel_PrepareAndSendData (uint32_t channel, array< int32_t >^ amplitude, array< uint64_t >^ duration, STG_DestinationEnumNet destType, bool doClear)
- StimulusDeviceDataAndUnrolledData ^ PrepareData (int channel, array< int32_t >^ amplitude, array< uint64 t >^ duration, STG DestinationEnumNet destType)
- void SendPreparedData (int channel, StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled, STG_DestinationEnumNet destType)
- SidebandData ^ CreateSideband (array< int32_t >^ StimulusActive, array< int32_t >^ Syncout, array< uint64 t >^ Duration, uint32 t Bit0Time, uint32 t Bit3Time, uint32 t Bit4Time)

Creates the Sideband Channel for the MEA2100 device.

void ClearMultiplexedData ()

Clears the Stimulation Memory in the STG device.

void SendMultiplexedData (array< uint16_t >^ data)

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

int GetMultiplexedDataChannelsInBlock ()

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData

• int GetDACResolution ()

Gets number of bits of the DAC resolution.

int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

int GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

int GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

• int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

void SetupTrigger (uint32_t first_trigger, array< uint32_t >^ channelmap, array< uint32_t >^ syncoutmap, array< uint32_t >^ repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

• void SetupTriggerSingle (uint32_t trigger, uint32_t channelmap, uint32_t syncoutmap, uint32_t repeat)

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

uint32_t GetTotalMemory ()

Get the total amount of memory available on the STG (all segments).

• uint32 t GetAvailableMemory ()

Get the amount of memory available in the currently selected segment of the STG.

• int GetNumberOfAnalogChannels ()

Get the number of STG channels.

Events

OnStgPollStatus^ PollStatusEvent

Additional Inherited Members

11.115.1 Constructor & Destructor Documentation

11.115.2 Member Function Documentation

Delete a Stimulus Pattern from STG memory

Parameters

```
channel specifies the channel to clear.
```

11.115.2.3 ClearMultiplexedData() void ClearMultiplexedData ()

Clears the Stimulation Memory in the STG device.

```
11.115.2.4 ClearSyncData() void ClearSyncData ( int channel )
```

Delete a Syncout Pattern from STG memory

channel	specifies the channel to clear.
---------	---------------------------------

Creates the Sideband Channel for the MEA2100 device.

Each datapoint is represented by an signed 32bit integer value. A value 0 means that the stimulation is active during that time. A value 1 means that the stimulation is not active during that time.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Parameters

SActive A list of datapoints which define weather the Stimulus is active or idle at that time as int32.

Parameters

Duration	A list of durations as uint64. The time is given in units of μs.
Bit0Time	Time in µs for which Bit 0 (Blanking) is to be extended.

Parameters

Bit3Time Time in μs for which Bit 3 (Stimulus Enable) is to be extende
--

Returns

Error Status. 0 on success.

11.115.2.6 ForceStatusEvent() void ForceStatusEvent ()

Force a status event. Force the DLL to create a PollMessage event and to call the pPollCallback function, even if no new status information is available.

11.115.2.7 GetAvailableMemory() uint32_t GetAvailableMemory ()

Get the amount of memory available in the currently selected segment of the STG.

Returns

The total memory available on the STG in bytes.

11.115.2.8 **GetCurrentRangeInNanoAmp()** int GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.115.2.9 GetCurrentResolutionInNanoAmp() int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

channel Channel which is	queried.
--------------------------	----------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.115.2.10 GetDACResolution() int GetDACResolution ()

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

11.115.2.11 GetMultiplexedDataChannelsInBlock() int GetMultiplexedDataChannelsInBlock ()

Gets the number of stimulus data channels to send per time slice. Might be greater than the number of configured channels. Fill unused channels with dummy data in SendMultiplexedData

$\textbf{11.115.2.12} \quad \textbf{GetNumberOfAnalogChannels()} \quad \texttt{int GetNumberOfAnalogChannels ()}$

Get the number of STG channels.

Returns

The number of STG channels.

11.115.2.13 GetTotalMemory() uint32_t GetTotalMemory ()

Get the total amount of memory available on the STG (all segments).

Returns

The total memory available on the STG in bytes.

11.115.2.14 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

channel	Channel which is queried.
---------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

11.115.2.15 GetVoltageResolutionInMicroVolt() int GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

Returns

The Voltage Resolution of the specified channel in Microvolts.

 ${\tt STG_DestinationEnumNet} \ \ \textit{destType} \ \)$ Prepare and append data to a given channel on the STG.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 μ s.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

channel The channel number to send dat	a to.
--	-------

amplitude	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
-----------	--

Parameters

duration	A list of durations in units of μs .
destType	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

Prepare and send data to a given channel on the STG. Previous data sent to that channel is erased first.

Each datapoint is represented by an signed 32bit integer value. When using voltage stimulation, the values are in multiple of 1 uV, thus the possible range is += 2000 V. When using current stimulation, the values are in multiple of 1 nA, this the possible range is += 2000 mA.

The duration is given as a list of 64 bit integers. Durations are given in units of μ s. The STG has a resolution of 20 us.

Blocks of data which should repeat can be defined by prepending such a block with an entry in the arrays where both amplitude and duration is zero. The end of such an block is marked by an entry where the duration is set to zero and the amplitude beeing set to the number of times the block should run. Blocks can be nested.

Parameters

channel	The channel number to send data to.

amplituda	A list of amplitudes in units of μV and nA in voltage and current mode, respectively.
ampiiluae	A list of amplitudes in units of Liv and ha in voltage and current mode, respectively.
,	, , ,

duration	A list of durations in units of μs .
destType	specifies wheather the data is for syncout, current or voltage stimulation.

Returns

Error Status. 0 on success.

```
11.115.2.18 PrepareData() StimulusDeviceDataAndUnrolledData ^ PrepareData (
    int channel,
    array< int32_t >^ amplitude,
    array< uint64_t >^ duration,
    STG_DestinationEnumNet destType )
```

```
11.115.2.19 SendMultiplexedData() void SendMultiplexedData ( array < uint16_t >^{\wedge} data )
```

Sends stimulus data in multiplexed form. All 16 bits words for the enabled DAC and digital channels are muxed together per time slice.

Parameters

```
data Array of data to be sent.
```

```
11.115.2.21 SendStart() void SendStart ( uint32_t triggermap )
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

triggermap A bitmap of triggers which will be started.
--

```
11.115.2.22 SendStop() [1/2] void SendStop ( uint32_t triggermap )
```

Stop some or all triggers of the STG.

Parameters

tri	ggermap	A bitmap of triggers which will be stopped.
-----	---------	---

```
11.115.2.23 SendStop() [2/2] void SendStop ( uint32_t triggermap, int options )
```

Stop some or all triggers of the STG.

Parameters

triggermap	A bitmap of triggers which will be stopped.
options	bitmap of options, currently only STOP_OPTION_SAVESTOP (0x80) is defined, which bypasses
	the stop commands when a syncout assossiated with a given sync-out has bit 1 (0x02) set. Can
	be used e.g. to prevent a stop while a biphasic stimulation pulse is active

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

first_trigger	The number of the first trigger to change.
---------------	--

channelmap	For each trigger, a bitmap of channels that belong to this trigger.

Parameters

syncoutmap	For each trigger, a bitmap of syncouts that belong to this trigger.
repeat	For each trigger, define the number of times this trigger should be repeated.

11.115.2.25 SetupTriggerSingle() void SetupTriggerSingle (

```
uint32_t trigger,
uint32_t channelmap,
uint32_t syncoutmap,
uint32_t repeat )
```

Configures the trigger settings for the STG. Note that all memory segments have their own trigger setting.

Parameters

trigger	The trigger to change.
---------	------------------------

Parameters

-6	A lating an of also and also the class at the state and
cnanneimap	A bitmap of channels that belong to this trigger.

Parameters

syncoutmap	A bitmap of syncouts that belong to this trigger.
repeat	The number of times this trigger should be repeated.

11.115.2.26 StartPoll() void StartPoll ()

Starts the interrupt fetching thread and delivers events

11.115.2.27 StopPoll() void StopPoll ()

Stops the interrupt fetching thread and delivers events


11.115.3 Event Documentation

11.115.3.1 PollStatusEvent OnStgPollStatus^ PollStatusEvent

11.116 CSw2to64DeviceNet Class Reference

The class to control the MCS-USB-Sw2to64 device.

Inheritance diagram for CSw2to64DeviceNet:

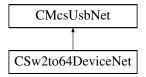

Public Member Functions

- CSw2to64DeviceNet ()
- ~CSw2to64DeviceNet ()
- unsigned short GetNumber ()

Gets the number of channels that can be switched in this box.

array< unsigned char > ^ GetChannels ()

Gets the current switch positions as char array.

void SetChannels (array< unsigned char >^ pattern)

Sets the switch positions from a char array.

unsigned char GetChannel (unsigned short index)

Gets one current switch position.

void SetChannel (unsigned short index, unsigned char pattern)

Sets one switch position.

Additional Inherited Members

11.116.1 Detailed Description

The class to control the MCS-USB-Sw2to64 device.

This class controls the settings of the MCS-USB-Sw2to64. The box has two inputs for signals. Each of the 64 outputs can be connected to one of the input signals, could be held open or connected ground. Valid switch states are 0, 1, 2 or 3 for each of the settings.

11.116.2 Constructor & Destructor Documentation

```
11.116.2.1 CSw2to64DeviceNet() CSw2to64DeviceNet ()
```

11.116.2.2 ~CSw2to64DeviceNet() ~CSw2to64DeviceNet ()

11.116.3 Member Function Documentation

Gets one current switch position.

Parameters

in	index	number of channel to read the switch position from
----	-------	--

Returns

switch position of desired channel

11.116.3.2 GetChannels() array<unsigned char> ^ GetChannels ()

Gets the current switch positions as char array.

Returns

array of char with the size of the number of channels, each char has the setting of a channel

11.116.3.3 GetNumber() unsigned short GetNumber ()

Gets the number of channels that can be switched in this box.

The box can have a different number of channels it can switch. Up to now usually 64 channels are returned

```
11.116.3.4 SetChannel() void SetChannel (
    unsigned short index,
    unsigned char pattern)
```

Sets one switch position.

in	index	number of channel to write the switch position to
in	pattern	switch position of the channel

11.116.3.5 SetChannels() void SetChannels ($array < unsigned char >^{\wedge} pattern$)

Sets the switch positions from a char array.


Parameters

in	pattern	array of char with the size of the number of channels, each char has the setting of a channel
----	---------	---

11.117 CTcxDeviceNet Class Reference

Class to control a Temperature Controller (TCX)

Inheritance diagram for CTcxDeviceNet:

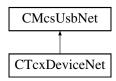

Public Member Functions

CTcxDeviceNet ()

Initializes a new instance of CTcxDeviceNet class.

- ∼CTcxDeviceNet ()
- unsigned int GetNumControlChannels ()

Gets the number of channels the device can control/regulate.

• unsigned int GetNumMeasureChannels ()

Gets the number of channels the device can measure.

• int GetValue (unsigned int channel)

Gets the temperate of the specified channel in units of 0.1 $^{\circ}$ C.

int GetValueHires (unsigned int channel)

Gets the temperate of the specified channel in units of 0.01 °C.

int GetHeaterTemp (unsigned int channel)

Gets the temperate of the specified heater in units of 0.1 °C.

• int GetHeaterLimit (unsigned int device)

Gets the temperate limit of the specified heater in units of 0.1 $^{\circ}\!\text{C}.$

• double GetMaxHeaterPowerMultiwell ()

queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

void SetMaxHeaterPowerMultiwell (double MaxPowerWatt)

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

bool GetHasThermocouple ()

Gets weather the device supports a thermocouple.

- bool GetEnableHeaterLimit (unsigned int device)
- bool GetEnableThermocouple (unsigned int device)
- TcxSensorTypeEnumNet GetSensorType (unsigned int device)
- String \(^\) GetUnit (unsigned int channel)
- unsigned int GetBoardTemp ()

Gets the temperate of the mainboard in units of 0.1 ℃.

- · unsigned int GetVolti (unsigned int channel)
- unsigned int GetNumDevices ()
- · void SetSetpoint (unsigned int channel, int sp)

Sets the target temperate of specified channel in units of 0.1 ℃.

- void SetDevice (unsigned int channel, int device)
- void SetOnOff (unsigned int channel, bool on)

Switches the specified channel on or off.

- void SetCalibration (unsigned int channel, int calib)
- void SetP (unsigned int device, int p_coeff)

Sets the P-coefficient of the specified device.

void SetI (unsigned int device, int i_coeff)

Sets the I-coefficient of the specified device.

void SetD (unsigned int device, int d_coeff)

Sets the D-coefficient of the specified device.

void SetMaxP (unsigned int device, int maxp)

Sets the maximum heater power of the specified device.

- void SetHeaterLimit (unsigned int device, int heater_limit)
- void SetEnableHeaterLimit (unsigned int device, bool enable)
- void SetEnableThermocouple (unsigned int device, bool enable)
- void SetSensorType (unsigned int device, TcxSensorTypeEnumNet type)
- void SetDevname (unsigned int device, String[^] Devicename)
- int GetSetpoint (unsigned int channel)

Gets the target temperate of specified channel in units of 0.1 °C.

- int GetDevice (unsigned int channel)
- int GetOnOff (unsigned int channel)

Gets if the specified channel is on or off.

- int GetCalibration (unsigned int channel)
- int GetP (unsigned int device)

Gets the P-coefficient of the specified device.

• int GetI (unsigned int device)

Gets the I-coefficient of the specified device.

int GetD (unsigned int device)

Gets the D-coefficient of the specified device.

• int GetMaxP (unsigned int device)

Gets the maximum heater power of the specified device.

- String \(^\) GetDevname (unsigned int device)
- TcxDeviceTypeEnumNet GetDeviceType ()
- int GetSetpointMin (unsigned int channel)
- int GetCalibrationMin (unsigned int channel)
- int GetPMin (unsigned int device)
- int GetIMin (unsigned int device)
- int GetDMin (unsigned int device)

- int GetMaxpMin (unsigned int device)
- int GetSetpointMax (unsigned int channel)
- int GetCalibrationMax (unsigned int channel)
- int GetPMax (unsigned int device)
- int GetIMax (unsigned int device)
- int GetDMax (unsigned int device)
- int GetMaxpMax (unsigned int device)
- int GetSetpointDecp (unsigned int channel)
- int GetCalibrationDecp (unsigned int channel)
- int GetPDecp (unsigned int device)
- int GetIDecp (unsigned int device)
- int GetDDecp (unsigned int device)
- int GetMaxpDecp (unsigned int device)
- int GetResX (unsigned int channel)
- int GetResS (unsigned int channel)
- int GetRes1 (unsigned int channel)
- int GetRes2 (unsigned int channel)
- int GetPwrSet (unsigned int channel)
- int GetPwrOut (unsigned int channel)
- int GetDuty (unsigned int channel)

Gets the duty cycle of the heating element.

int GetUOut (unsigned int channel)

Gets the voltage on the heating element.

• int GetlOut (unsigned int channel)

Gets the current through the heating element.

• int GetROut (unsigned int channel)

Gets the resistance of the heating element.

int GetPOut (unsigned int channel)

Gets the output power of the heating element.

- · int GetCurrent (unsigned int channel)
- int GetThermocoupleTemp (unsigned int channel)
- int GetThermocoupleTempAbs (unsigned int channel)
- int GetThermocoupleReferenceTemp (unsigned int channel)
- unsigned int GetThermocoupleNanovoltPerKelvin (unsigned int channel)

Gets the proportional constant for the thermocouple.

void SetThermocoupleNanovoltPerKelvin (unsigned int channel, unsigned int value)

Sets the proportional constant for the thermocouple.

- int GetThermocoupleCalibration (unsigned int channel)
- void CalibrateThermocouple (unsigned int channel)
- void SetDeviceType (TcxDeviceTypeEnumNet devicetype)
- void FactoryReset ()

Additional Inherited Members

11.117.1 Detailed Description

Class to control a Temperature Controller (TCX)

11.117.2 Constructor & Destructor Documentation

```
11.117.2.1 CTcxDeviceNet() CTcxDeviceNet ()
Initializes a new instance of CTcxDeviceNet class.
11.117.2.2 ~CTcxDeviceNet() ~CTcxDeviceNet ()
11.117.3 Member Function Documentation
11.117.3.1 CalibrateThermocouple() void CalibrateThermocouple (
             unsigned int channel )
11.117.3.2 FactoryReset() void FactoryReset ( )
11.117.3.3 GetBoardTemp() unsigned int GetBoardTemp ( )
Gets the temperate of the mainboard in units of 0.1 °C.
11.117.3.4 GetCalibration() int GetCalibration (
             unsigned int channel )
11.117.3.5 GetCalibrationDecp() int GetCalibrationDecp (
             unsigned int channel )
11.117.3.6 GetCalibrationMax() int GetCalibrationMax (
             unsigned int channel )
11.117.3.7 GetCalibrationMin() int GetCalibrationMin (
             unsigned int channel )
```

```
11.117.3.8 GetCurrent() int GetCurrent (
             unsigned int channel )
11.117.3.9 GetD() int GetD (
             unsigned int device )
Gets the D-coefficient of the specified device.
11.117.3.10 GetDDecp() int GetDDecp (
             unsigned int device )
11.117.3.11 GetDevice() int GetDevice (
             unsigned int channel )
11.117.3.12 GetDeviceType() TcxDeviceTypeEnumNet GetDeviceType ( )
11.117.3.13 GetDevname() String ^{\wedge} GetDevname (
             unsigned int device )
11.117.3.14 GetDMax() int GetDMax (
             unsigned int device )
11.117.3.15 GetDMin() int GetDMin (
             unsigned int device )
11.117.3.16 GetDuty() int GetDuty (
             unsigned int channel )
```

Gets the duty cycle of the heating element.

Generated by Doxygen

Returns

The duty cycle in percent, the value of 320 * 64 corresponds to 100 %.

```
11.117.3.17 GetEnableHeaterLimit() bool GetEnableHeaterLimit ( unsigned int device )
```

```
11.117.3.18 GetEnableThermocouple() bool GetEnableThermocouple ( unsigned int device )
```

11.117.3.19 GetHasThermocouple() bool GetHasThermocouple ()

Gets weather the device supports a thermocouple.

Gets the temperate limit of the specified heater in units of 0.1 ℃.

Gets the temperate of the specified heater in units of 0.1 °C.

```
11.117.3.22 Getl() int GetI ( unsigned int device )
```

Gets the I-coefficient of the specified device.

```
11.117.3.23 GetIDecp() int GetIDecp (
              unsigned int device )
11.117.3.24 GetIMax() int GetIMax (
             unsigned int device )
11.117.3.25 GetlMin() int GetIMin (
             unsigned int device )
11.117.3.26 GetIOut() int GetIOut (
              unsigned int channel )
Gets the current through the heating element.
Parameters
 channel The channel number.
Returns
     The current in units of mA.
11.117.3.27 GetMaxHeaterPowerMultiwell() double GetMaxHeaterPowerMultiwell ( )
queries the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W
11.117.3.28 GetMaxP() int GetMaxP (
              unsigned int device )
Gets the maximum heater power of the specified device.
```

11.117.3.29 **GetMaxpDecp()** int GetMaxpDecp (unsigned int *device*)

```
11.117.3.30 GetMaxpMax() int GetMaxpMax (
             unsigned int device )
11.117.3.31 GetMaxpMin() int GetMaxpMin (
             unsigned int device )
11.117.3.32 GetNumControlChannels() unsigned int GetNumControlChannels ( )
Gets the number of channels the device can control/regulate.
11.117.3.33 GetNumDevices() unsigned int GetNumDevices ()
11.117.3.34 GetNumMeasureChannels() unsigned int GetNumMeasureChannels ( )
Gets the number of channels the device can measure.
11.117.3.35 GetOnOff() int GetOnOff (
             unsigned int channel )
Gets if the specified channel is on or off.
11.117.3.36 GetP() int GetP (
             unsigned int device )
Gets the P-coefficient of the specified device.
11.117.3.37 GetPDecp() int GetPDecp (
             unsigned int device )
11.117.3.38 GetPMax() int GetPMax (
             unsigned int device )
11.117.3.39 GetPMin() int GetPMin (
             unsigned int device )
11.117.3.40 GetPOut() int GetPOut (
             unsigned int channel )
Gets the output power of the heating element.
```

channel The channel number	er.
----------------------------	-----

Returns

The resistance in units of mW.

```
\textbf{11.117.3.41} \quad \textbf{GetPwrOut()} \quad \texttt{int GetPwrOut} \ \ (
              unsigned int channel )
11.117.3.42 GetPwrSet() int GetPwrSet (
              unsigned int channel )
11.117.3.43 GetRes1() int GetRes1 (
              unsigned int channel )
11.117.3.44 GetRes2() int GetRes2 (
              unsigned int channel )
11.117.3.45 GetResS() int GetResS (
              unsigned int channel )
11.117.3.46 GetResX() int GetResX (
             unsigned int channel )
11.117.3.47 GetROut() int GetROut (
              unsigned int channel )
```

Gets the resistance of the heating element.

channel The channel number.	
-----------------------------	--

Returns

The resistance in units of 0.1 Ohm.

```
11.117.3.48 GetSensorType() TcxSensorTypeEnumNet GetSensorType ( unsigned int device )
```

```
11.117.3.49 GetSetpoint() int GetSetpoint (
          unsigned int channel )
```

Gets the target temperate of specified channel in units of 0.1 ℃.

```
11.117.3.53 GetThermocoupleCalibration() int GetThermocoupleCalibration ( unsigned int channel)
```

```
11.117.3.54 GetThermocoupleNanovoltPerKelvin() unsigned int GetThermocoupleNanovoltPerKelvin ( unsigned int channel)
```

Gets the proportional constant for the thermocouple.

channel	Thermocouple channel number.
---------	------------------------------

Returns

The proportional constant in Nanovolt per Kelvin.

```
11.117.3.58 GetUnit() String ^ GetUnit (
unsigned int channel)
```

```
11.117.3.59 GetUOut() int GetUOut ( unsigned int channel)
```

Gets the voltage on the heating element.

Parameters

channel	The channel number.
or iai ii ioi	The charmer hamben

Returns

The voltage in units of mV.

```
11.117.3.60 GetValue() int GetValue (
unsigned int channel)
```

Gets the temperate of the specified channel in units of 0.1 ℃.

```
11.117.3.61 GetValueHires() int GetValueHires ( unsigned int channel)
```

Gets the temperate of the specified channel in units of 0.01 $^{\circ}\!\text{C}.$

```
11.117.3.62 GetVolti() unsigned int GetVolti ( unsigned int channel )
```

```
11.117.3.63 SetCalibration() void SetCalibration (
    unsigned int channel,
    int calib )
```

```
11.117.3.64 SetD() void SetD (
unsigned int device,
int d_coeff)
```

Sets the D-coefficient of the specified device.

```
11.117.3.65 SetDevice() void SetDevice (
unsigned int channel,
int device )
```

```
11.117.3.67 SetDevname() void SetDevname (
unsigned int device,
String^ Devicename)
```

```
11.117.3.72 SetMaxHeaterPowerMultiwell() void SetMaxHeaterPowerMultiwell (

double MaxPowerWatt )
```

sets the max. heater power that the Multiwell temperature controller will apply; unit: W; useful range: 5.2W..7.6W

```
11.117.3.73 SetMaxP() void SetMaxP (
unsigned int device,
int maxp )
```

Sets the maximum heater power of the specified device.

```
11.117.3.74 SetOnOff() void SetOnOff (
     unsigned int channel,
     bool on )
```

Switches the specified channel on or off.

Parameters

channel The channel number.

```
11.117.3.75 SetP() void SetP ( unsigned int device, int p\_coeff )
```

Sets the P-coefficient of the specified device.

```
11.117.3.77 SetSetpoint() void SetSetpoint ( unsigned int channel, int sp )
```

Sets the target temperate of specified channel in units of 0.1 ℃.

```
11.117.3.78 SetThermocoupleNanovoltPerKelvin() void SetThermocoupleNanovoltPerKelvin (
    unsigned int channel,
    unsigned int value )
```

Sets the proportional constant for the thermocouple.


Parameters

channel	Thermocouple channel number.
value	Proportinal constant in Nanovolt per Kelvin.

11.118 CTEERFunctionNet Class Reference

CTEERFunctionNet is the class to control the TEER device

Inheritance diagram for CTEERFunctionNet:

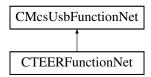

Public Member Functions

 uint32_t GetRotaryPositionCode () gets the rotary position code

 CTEERFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pTEERFunctionPointer← Container) Initializes a new instance of the CTEERFunctionNet class. CTEERFunctionNet (CMcsUsbNet[^] mcsusb) virtual ~CTEERFunctionNet () • !CTEERFunctionNet () • uint32 t GetPeriod us () gets the period of TEER stimulation in us void SetPeriod_us (uint32_t period_us) sets the period of TEER stimulation in us uint32 t GetAmplitude nA () gets TEER stimulation amplitude in nA void SetAmplitude_nA (uint32_t Amplitude_nA) sets TEER stimulation amplitude in nA TeerWaveformEnumNet GetWaveform () gets TEER stimulation waveform (sine/rect) void SetWaveform (TeerWaveformEnumNet Waveform) sets TEER stimulation waveform (sine/rect) TeerClampModeEnumNet GetClampMode () gets TEER clamp mode (voltage/current) void SetClampMode (TeerClampModeEnumNet ClampMode) sets TEER clamp mode (voltage/current) void StartSampling (uint32_t NumberOfCycles) starts TEER stimulation (duration: n cycles) and samples during last cycle · void StopSampling () stops TEER stimulation and sampling uint32 t IsSamplingFinished () returns false iff stimulation/sampling is going on, otherwise true void SetControllerParams (uint32 t P, uint32 t I, uint32 t D) sets PID controller parameters for voltage clamp mode void GetControllerParams ([System::Runtime::InteropServices::Out]uint32 t% P, [System::Runtime::← InteropServices::Out]uint32_t% I, [System::Runtime::InteropServices::Out]uint32_t% D) gets PID controller parameters for voltage clamp mode array< int32 t > ^ GetSampleBufferChunk (int Buffer Length) private function to query max. 100 bytes of sample buffer; called internally array< int32_t > ^ GetSampleVoltageBuffer_uV (int Buffer_Length) returns voltage sample buffer (max. 500 values); unit: uV uint32 t GetMaxChunkSize Byte () private function to be called internally only uint32_t GetBytesPerSample () private function to be called internally only uint32 t GetNumberOfAvailableSamples () private function to be called internally only void SetBufferIndex (uint32 t NewBufferIndex) pre-selects sample buffer to be tranferred by GetSampleVoltageBuffer_uV() uint32_t GetAdapterCode () gets the adapter code

```
    void SetExternalLED (uint32_t NewState)

      sets the external LED

    void SetCurrentEnable (bool NewCurrentEnable)

      when disabled, no current will flow through chamber

    bool GetCurrentEnable ()

      when disabled, no current will flow through chamber
• int32 t GetUptimeSeconds ()
      returns time in seconds since device was powered up

    void StartInternalCalibration ()

      starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call
• bool IsInternalCalibrationFinished ()
      queries whether internal calibration has finished

    int GetDacZero ()

      returns DAC-offset (result of internal calibration); use to check for plausibility only

    void CancelInternalCalibration ()

      in case the internal calibration "hangs", this will cancel it
• void SetLiquidResistance (int32_t NewLiquidResistance_Ohm)
      sets the resistance of the liquid in ohms
• int32_t GetLiquidResistance ()
      gets the resitance of the liquid in ohms

    int GetScaleFactorU1 ()

      returns U1 scale factor times 10<sup>^</sup>6 (result of internal calibration)

    int GetScaleFactorU2 ()

      returns U2 scale factor times 10<sup>6</sup> (result of internal calibration)

    int GetAdcOffsetU1 ()

      returns ADC offset of U1 channel (result of internal calibration)
• int GetAdcOffsetU2 ()
      returns ADC offset of U2 channel (result of internal calibration)
```

returns number of times (since bootup) sample memory got overwritten

Additional Inherited Members

int GetSampleRate ()

int GetFrameErrorCounter ()

returns sample rate in Hz

11.118.1 Detailed Description

CTEERFunctionNet is the class to control the TEER device

11.118.2 Constructor & Destructor Documentation

```
11.118.2.1 CTEERFunctionNet() [1/2] CTEERFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pTEERFunctionPointerContainer)
```

Initializes a new instance of the CTEERFunctionNet class.

```
11.118.2.2 CTEERFunctionNet() [2/2] CTEERFunctionNet (
              CMcsUsbNet^ mcsusb )
11.118.2.3 ~CTEERFunctionNet() virtual ~CTEERFunctionNet ( ) [virtual]
11.118.2.4 "!CTEERFunctionNet() !CTEERFunctionNet ( )
11.118.3 Member Function Documentation
11.118.3.1 CancelInternalCalibration() void CancelInternalCalibration ( )
in case the internal calibration "hangs", this will cancel it
11.118.3.2 GetAdapterCode() uint32_t GetAdapterCode ( )
gets the adapter code
Returns
     the adapter code
11.118.3.3 GetAdcOffsetU1() int GetAdcOffsetU1 ( )
returns ADC offset of U1 channel (result of internal calibration)
Returns
     the ADC offset for U1
11.118.3.4 GetAdcOffsetU2() int GetAdcOffsetU2 ( )
returns ADC offset of U2 channel (result of internal calibration)
Returns
     the ADC offset for U2
```

```
11.118.3.5 GetAmplitude_nA() uint32_t GetAmplitude_nA ( )
```

gets TEER stimulation amplitude in nA

Returns

current stimulation amplitude in nA

```
11.118.3.6 GetBytesPerSample() uint32_t GetBytesPerSample ( )
```

private function to be called internally only

Returns

not documented

```
11.118.3.7 GetClampMode() TeerClampModeEnumNet GetClampMode ( )
```

gets TEER clamp mode (voltage/current)

Returns

current TEER clamp mode

```
11.118.3.8 GetControllerParams() void GetControllerParams (
```

```
[System::Runtime::InteropServices::Out] uint32_t% P,
[System::Runtime::InteropServices::Out] uint32_t% I,
[System::Runtime::InteropServices::Out] uint32_t% D)
```

gets PID controller parameters for voltage clamp mode

Parameters

Р	the P value
1	the I value
D	the D value

11.118.3.9 GetCurrentEnable() bool GetCurrentEnable ()

when disabled, no current will flow through chamber

Returns

not documented

false when disabled, true when enabled

```
11.118.3.10 GetDacZero() int GetDacZero ( )
returns DAC-offset (result of internal calibration); use to check for plausibility only
Returns
     the DAC offset
11.118.3.11 GetFrameErrorCounter() int GetFrameErrorCounter ( )
returns number of times (since bootup) sample memory got overwritten
Returns
     the number of errors
11.118.3.12 GetLiquidResistance() int32_t GetLiquidResistance ( )
gets the resitance of the liquid in ohms
Returns
     the resistance in ohms
11.118.3.13 GetMaxChunkSize_Byte() uint32_t GetMaxChunkSize_Byte ( )
private function to be called internally only
Returns
```

```
11.118.3.14 GetNumberOfAvailableSamples() uint32_t GetNumberOfAvailableSamples ( )
private function to be called internally only
Returns
     not documented
11.118.3.15 GetPeriod_us() uint32_t GetPeriod_us ( )
gets the period of TEER stimulation in us
Returns
     the period in us
11.118.3.16 GetRotaryPositionCode() uint32_t GetRotaryPositionCode ( )
gets the rotary position code
Returns
     the rotary position code
11.118.3.17 GetSampleBufferChunk() array<int32_t> ^{\land} GetSampleBufferChunk (
              int Buffer_Length )
private function to query max. 100 bytes of sample buffer; called internally
Parameters
 Buffer_Length
                 The maximal length of Buffer.
Returns
     not documented
```

Generated by Doxygen

returns sample rate in Hz

11.118.3.18 GetSampleRate() int GetSampleRate ()

_	-		
О	A *1	1114	20
п	eu	ui	115

the sample rate in Hz

```
11.118.3.19 GetSampleVoltageBuffer_uV() array<int32_t> ^{\land} GetSampleVoltageBuffer_uV ( int Buffer_Length )
```

returns voltage sample buffer (max. 500 values); unit: uV

Parameters

Buffer_Length	The maximal length of Buffer.
---------------	-------------------------------

Returns

the voltage sample buffer

11.118.3.20 GetScaleFactorU1() int GetScaleFactorU1 ()

returns U1 scale factor times 10⁶ (result of internal calibration)

Returns

the U1 scale factor

$\textbf{11.118.3.21} \quad \textbf{GetScaleFactorU2()} \quad \texttt{int GetScaleFactorU2 ()} \\$

returns U2 scale factor times $10^{\circ}6$ (result of internal calibration)

Returns

the U2 scale factor

11.118.3.22 GetUptimeSeconds() int32_t GetUptimeSeconds ()

returns time in seconds since device was powered up

Returns

seconds since power-on

```
11.118.3.23 GetWaveform() TeerWaveformEnumNet GetWaveform ( )
gets TEER stimulation waveform (sine/rect)
Returns
     waveform enum
11.118.3.24 IsInternalCalibrationFinished() bool IsInternalCalibrationFinished ( )
queries whether internal calibration has finished
Returns
     true if calibration has finished
11.118.3.25 IsSamplingFinished() uint32_t IsSamplingFinished ( )
returns false iff stimulation/sampling is going on, otherwise true
Returns
     true if sampling is finished
11.118.3.26 SetAmplitude_nA() void SetAmplitude_nA (
              uint32_t Amplitude_nA )
sets TEER stimulation amplitude in nA
Parameters
 Amplitude_nA
                 new stimulation amplitude in nA
```

```
11.118.3.27 SetBufferIndex() void SetBufferIndex ( uint32_t NewBufferIndex )
```

pre-selects sample buffer to be tranferred by GetSampleVoltageBuffer_uV()

Parameters

NewBufferIndex	0 - chamber voltage; 1 - compliance voltage

sets TEER clamp mode (voltage/current)

Parameters

sets PID controller parameters for voltage clamp mode

Parameters

Р	the P value
1	the I value
D	the D value

11.118.3.30 SetCurrentEnable() void SetCurrentEnable (bool NewCurrentEnable)

when disabled, no current will flow through chamber

Parameters

NewCurrentEnable disabled when false, enabled when true

```
11.118.3.31 SetExternalLED() void SetExternalLED ( uint32_t NewState )
```

sets the external LED

Parameters

NewState state

```
11.118.3.32 SetLiquidResistance() void SetLiquidResistance ( int32_t NewLiquidResistance_Ohm )
```

sets the resistance of the liquid in ohms

Parameters

NewLiquidResistance_Ohm	the resistance in ohms
-------------------------	------------------------

```
11.118.3.33 SetPeriod_us() void SetPeriod_us ( uint32_t period_us )
```

sets the period of TEER stimulation in us

Parameters

```
11.118.3.34 SetWaveform() void SetWaveform (

TeerWaveformEnumNet Waveform )
```

sets TEER stimulation waveform (sine/rect)

Parameters

Waveform | waveform enum

11.118.3.35 StartInternalCalibration() void StartInternalCalibration ()

starts determination of internal DAC-offset; result is used internally; NON-BLOCKING call

```
11.118.3.36 StartSampling() void StartSampling ( uint32_t NumberOfCycles )
```

starts TEER stimulation (duration: n cycles) and samples during last cycle

Parameters


NumberOfCvcles	number of cycles (sine or rect) to output (0 - loop forever)

11.118.3.37 StopSampling() void StopSampling ()

stops TEER stimulation and sampling

11.119 CTEERMachineDeviceNet Class Reference

Inheritance diagram for CTEERMachineDeviceNet:

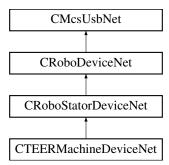

Public Member Functions

- CTEERMachineDeviceNet ()
- ∼CTEERMachineDeviceNet ()

Properties

• CTEERFunctionNet [get]

Additional Inherited Members

11.119.1 Constructor & Destructor Documentation

11.119.1.1 CTEERMachineDeviceNet() CTEERMachineDeviceNet ()

11.119.1.2 ~CTEERMachineDeviceNet() ~CTEERMachineDeviceNet ()


11.119.2 Property Documentation

11.119.2.1 TEERFunctionNet CTEERFunctionNet [get]

11.120 CUsbDeviceConfigurationFunctionNet Class Reference

CUsbDeviceConfigurationFunctionNet is the class to configure the USB firmware

Inheritance diagram for CUsbDeviceConfigurationFunctionNet:

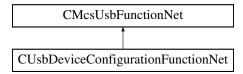

Public Member Functions

- CUsbDeviceConfigurationFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] p
 UsbDeviceConfigurationFunctionPointerContainer)
 - Initializes a new instance of the CUsbDeviceConfigurationFunctionNet class.
- CUsbDeviceConfigurationFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual \sim CUsbDeviceConfigurationFunctionNet ()
- !CUsbDeviceConfigurationFunctionNet ()
- void SetDeviceName (String^ name)

sets the USB device name for configurable devices

void SetDeviceId (ProductIdEnumNet id)

sets the USB device name for configurable devices

Additional Inherited Members

11.120.1 Detailed Description

CUsbDeviceConfigurationFunctionNet is the class to configure the USB firmware

11.120.2 Constructor & Destructor Documentation

```
11.120.2.1 CUsbDeviceConfigurationFunctionNet() [1/2] CUsbDeviceConfigurationFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pUsbDeviceConfigurationFunctionPointerContainer)
```


Initializes a new instance of the CUsbDeviceConfigurationFunctionNet class.

```
11.120.2.2 CUsbDeviceConfigurationFunctionNet() [2/2] CUsbDeviceConfigurationFunctionNet (
                                                                                                                        CMcsUsbNet^ mcsusb )
\textbf{11.120.2.3} \quad \sim \textbf{CUsbDeviceConfigurationFunctionNet()} \quad \text{virtual } \sim \texttt{CUsbDeviceConfigurationFunctionNet} \quad \textbf{(} 
 ) [virtual]
11.120.2.4 "!CUsbDeviceConfigurationFunctionNet() !CUsbDeviceConfigurationFunctionNet ( )
 11.120.3 Member Function Documentation
11.120.3.1 SetDeviceId() void SetDeviceId (
                                                                                                                       ProductIdEnumNet id )
 sets the USB device name for configurable devices
  Parameters
              id
11.120.3.2 SetDeviceName() void SetDeviceName (
                                                                                                                       String^{\wedge} name )
sets the USB device name for configurable devices
  Parameters
            name
```

11.121 CUsbExceptionNet Class Reference

Exception class that is thrown in case of an USB error.

Inheritance diagram for CUsbExceptionNet:

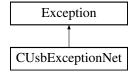

Public Member Functions

CUsbExceptionNet (uint32_t status)

Constructor of a CUsbException.

• CUsbExceptionNet (uint32_t status, String^ message)

Properties

```
• uint32_t Status [get]
```

11.121.1 Detailed Description

Exception class that is thrown in case of an USB error.

11.121.2 Constructor & Destructor Documentation

```
11.121.2.1 CUsbExceptionNet() [1/2] CUsbExceptionNet ( uint32_t status )
```

Constructor of a CUsbException.

Parameters

status the status number

11.121.3 Property Documentation

```
11.121.3.1 Status uint32_t Status [get]
```

11.122 CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet Class Reference

Public Member Functions

CVoltageRangeInfoNet (int vr, String[^] vrString)

Public Attributes

- int VoltageRangeInMicroVolt
- String ^ VoltageRangeDisplayStringMilliVolt

11.122.1 Constructor & Destructor Documentation

```
11.122.1.1 CVoltageRangeInfoNet() CVoltageRangeInfoNet ( int vr, String^{\wedge} vrString)
```


11.122.2 Member Data Documentation

11.122.2.1 VoltageRangeDisplayStringMilliVolt String ^ VoltageRangeDisplayStringMilliVolt

 $\textbf{11.122.2.2} \quad \textbf{VoltageRangeInMicroVolt} \quad \texttt{int VoltageRangeInMicroVolt}$

11.123 CW2100_FunctionNet Class Reference

Inheritance diagram for CW2100_FunctionNet:

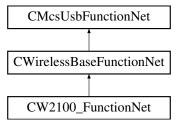

Classes

struct AudioChannelsNet

Public Member Functions

- CW2100_FunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] w2100_Function
 — PointerContainer)
- CW2100_FunctionNet (CMcsUsbNet[^] mcsusb)
- array< HeadStageIDType[^]> [^] GetAvailableHeadstages (unsigned int max_length)
- void SelectHeadstage (unsigned int IDorEntry, int TimeSlotNr)
- void DeselectHeadstage (int TimeSlotNr)
- void DeselectAllHeadstages ()
- HeadStageIDTypeState ^ GetSelectedHeadstageState (int TimeSlotNr)
- BatteryState ^ GetBatteryState (int TimeSlotNr)
- System::String ^ GetUserDefinedName (unsigned short ID)
- System::String ^ GetUserDefinedNameFromSelectedHS (int TimeSlotNr)
- System::String \(^\) GetUserDefinedNameCache (unsigned short ID)
- W2100_StimulusParametersNet ^ GetStiumlusParameters (unsigned short ID)
- W2100 StimulusParametersNet ^ GetStimulusParametersFromSelectedHS (int TimeSlotNr)
- W2100 StimulusParametersNet ^ GetStimulusParametersCache (unsigned int typeValue)
- uint32_t GetStimulusParametersCache (unsigned int typeValue, [System::Runtime::InteropServices::
 — Out]W2100_StimulusParametersNet^% StimulusParameters)
- void SetSelectedChannels (array< BYTE >^ channels, int TimeSlotNr)
- array< BYTE > ^ GetSelectedChannels (int TimeSlotNr)
- void SetMultiHeadstageMode (bool Mode)
- bool GetMultiHeadstageMode ()
- void SetHeadstageSamplingActive (bool Active, int TimeSlotNr)
- bool GetHeadstageSamplingActive (int TimeSlotNr)
- void SetHeadstageToSleep (unsigned int Sleep16ms, int TimeSlotNr)
- void SetHeadstageOnOff (unsigned short On, int TimeSlotNr)
- unsigned short GetHeadstageOnOff (int TimeSlotNr)
- unsigned int GetAnalogOutChannel ([System::Runtime::InteropServices::Out]int % automatic, unsigned short index)
- void SetAnalogOutChannel (int automatic, unsigned short index, unsigned int Channel)
- array< unsigned int > ^ GetAnalogOutFilter ([System::Runtime::InteropServices::Out]int % automatic)
- void SetAnalogOutFilter (int automatic, array< unsigned int >^ Coeffs)
- AnalogOut_DAC_Range_EnumNet GetDacRange ()
- void SetDacRange (AnalogOut_DAC_Range_EnumNet range)
- CFilterPropertyNet ^ GetFilterProperty (W2100DacqGroupChannelEnumNet GroupID, unsigned int index)
- array< CFilterPropertyNet[^]> [^] GetFilterProperties (W2100DacqGroupChannelEnumNet GroupID)
- void SetAccelGyroEnabled (W2100_Accel_Gyro_Select_EnumNet enable, int TimeSlotNr)
- W2100_Accel_Gyro_Select_EnumNet GetAccelGyroEnabled (int TimeSlotNr)
- void SetAccelGyroDesiredRate (int rate, int TimeSlotNr)
- int GetAccelGyroDesiredRate (int TimeSlotNr)
- int GetAccelGyroCurrentRate (int TimeSlotNr)
- void SetAccelRange (int range, int TimeSlotNr)
- int GetAccelRange (int TimeSlotNr)
- void SetGyroRange (int range, int TimeSlotNr)
- int GetGyroRange (int TimeSlotNr)
- void SetAudioChannels (array< AudioChannelsNet[^]>[^] channels)
- array< AudioChannelsNet[^]> [^] GetAudioChannels ()
- unsigned int GetPicFirmwareType (int TimeSlotNr)
- unsigned int GetFPGAFirmwareType (int TimeSlotNr)

Static Public Member Functions

- static void ClearUserDefinedNameCache ()
- static void ClearUserDefinedNameCache (unsigned short ID)
- static void ClearStimulusParametersCache ()
- static void ClearStimulusParametersCache (unsigned short ID)

Properties

- CW2100_StimulatorFunctionNet^ Stimulator [get]
- CPulseGeneratorFunctionNet^ PulseGenerator [get]

Additional Inherited Members

11.123.1 Constructor & Destructor Documentation

```
11.123.1.2 CW2100_FunctionNet() [2/2] CW2100_FunctionNet (
CMcsUsbNet^ mcsusb )
```

11.123.2 Member Function Documentation

```
11.123.2.1 ClearStimulusParametersCache() [1/2] static void ClearStimulusParametersCache ( ) [static]
```

```
11.123.2.2 ClearStimulusParametersCache() [2/2] static void ClearStimulusParametersCache ( unsigned short ID ) [static]
```

11.123.2.3 ClearUserDefinedNameCache() [1/2] static void ClearUserDefinedNameCache () [static]

```
11.123.2.4 ClearUserDefinedNameCache() [2/2] static void ClearUserDefinedNameCache (
              unsigned short ID ) [static]
11.123.2.5 DeselectAllHeadstages() void DeselectAllHeadstages ()
11.123.2.6 DeselectHeadstage() void DeselectHeadstage (
              int TimeSlotNr )
11.123.2.7 GetAccelGyroCurrentRate() int GetAccelGyroCurrentRate (
              int TimeSlotNr )
11.123.2.8 GetAccelGyroDesiredRate() int GetAccelGyroDesiredRate (
              int TimeSlotNr )
\textbf{11.123.2.9} \quad \textbf{GetAccelGyroEnabled()} \quad \texttt{W2100\_Accel\_Gyro\_Select\_EnumNet} \quad \texttt{GetAccelGyroEnabled} \quad \textbf{(}
              int TimeSlotNr )
11.123.2.10 GetAccelRange() int GetAccelRange (
              int TimeSlotNr )
\textbf{11.123.2.11} \quad \textbf{GetAnalogOutChannel()} \quad \texttt{unsigned int GetAnalogOutChannel ()}
              [System::Runtime::InteropServices::Out] int % automatic,
              unsigned short index)
11.123.2.12 GetAnalogOutFilter() array<unsigned int> ^ GetAnalogOutFilter (
              [System::Runtime::InteropServices::Out] int % automatic )
11.123.2.13 GetAudioChannels() array<AudioChannelsNet^> ^ GetAudioChannels ( )
```

```
11.123.2.14 GetAvailableHeadstages() array<HeadStageIDType^> ^ GetAvailableHeadstages (
              unsigned int max_length )
11.123.2.15 GetBatteryState() BatteryState ^ GetBatteryState (
              int TimeSlotNr )
11.123.2.16 GetDacRange() AnalogOut_DAC_Range_EnumNet GetDacRange ( )
11.123.2.17 GetFilterProperties() array<CFilterPropertyNet^> ^ GetFilterProperties (
              W2100DacqGroupChannelEnumNet GroupID )
11.123.2.18 GetFilterProperty() CFilterPropertyNet ^ GetFilterProperty (
              W2100DacqGroupChannelEnumNet GroupID,
              unsigned int index )
11.123.2.19 GetFPGAFirmwareType() unsigned int GetFPGAFirmwareType (
              int TimeSlotNr )
11.123.2.20 GetGyroRange() int GetGyroRange (
              int TimeSlotNr )
{\bf 11.123.2.21} \quad {\bf GetHeadstageOnOff()} \quad {\tt unsigned short GetHeadstageOnOff} \end{\cite{Approximation}} \label{thm:condition}
              int TimeSlotNr )
11.123.2.22 GetHeadstageSamplingActive() bool GetHeadstageSamplingActive (
              int TimeSlotNr )
\textbf{11.123.2.23} \quad \textbf{GetMultiHeadstageMode()} \quad \texttt{bool GetMultiHeadstageMode ()} \\
```

```
11.123.2.24 GetPicFirmwareType() unsigned int GetPicFirmwareType (
             int TimeSlotNr )
11.123.2.25 GetSelectedChannels() array<BYTE> ^ GetSelectedChannels (
             int TimeSlotNr )
11.123.2.26 GetSelectedHeadstageState() HeadStageIDTypeState ^ GetSelectedHeadstageState (
             int TimeSlotNr )
11.123.2.27 GetStimulusParametersCache() [1/2] W2100_StimulusParametersNet ^ GetStimulus↔
ParametersCache (
            unsigned int typeValue )
11.123.2.28 GetStimulusParametersCache() [2/2] uint32_t GetStimulusParametersCache (
            unsigned int typeValue,
             [System::Runtime::InteropServices::Out] W2100_StimulusParametersNet^% Stimulus↔
Parameters )
11.123.2.29 GetStimulusParametersFromSelectedHS() w2100_StimulusParametersNet ^ GetStimulus↔
ParametersFromSelectedHS (
            int TimeSlotNr )
11.123.2.30 GetStiumlusParameters() W2100_StimulusParametersNet ^ GetStiumlusParameters (
             unsigned short ID )
11.123.2.31 GetUserDefinedName() System::String ^ GetUserDefinedName (
             unsigned short ID )
11.123.2.32 GetUserDefinedNameCache() [1/2] System::String ^ GetUserDefinedNameCache (
             unsigned short ID )
```

```
11.123.2.33 GetUserDefinedNameCache() [2/2] uint32_t GetUserDefinedNameCache (
             unsigned short ID,
             [System::Runtime::InteropServices::Out] System::String^{\wedge}% Name )
11.123.2.34 GetUserDefinedNameFromSelectedHS() System::String ^ GetUserDefinedNameFrom←
SelectedHS (
             int TimeSlotNr )
11.123.2.35 SelectHeadstage() void SelectHeadstage (
             unsigned int IDorEntry,
             int TimeSlotNr )
11.123.2.36 SetAccelGyroDesiredRate() void SetAccelGyroDesiredRate (
             int rate,
             int TimeSlotNr )
11.123.2.37 SetAccelGyroEnabled() void SetAccelGyroEnabled (
             W2100_Accel_Gyro_Select_EnumNet enable,
             int TimeSlotNr )
11.123.2.38 SetAccelRange() void SetAccelRange (
             int range,
             int TimeSlotNr )
11.123.2.39 SetAnalogOutChannel() void SetAnalogOutChannel (
             int automatic,
             unsigned short index,
             unsigned int Channel )
11.123.2.40 SetAnalogOutFilter() void SetAnalogOutFilter (
             int automatic,
             array< unsigned int >^{\land} Coeffs )
```

```
11.123.2.41 SetAudioChannels() void SetAudioChannels (
            11.123.2.42 SetDacRange() void SetDacRange (
            AnalogOut_DAC_Range_EnumNet range )
11.123.2.43 SetGyroRange() void SetGyroRange (
            int range,
            int TimeSlotNr )
11.123.2.44 SetHeadstageOnOff() void SetHeadstageOnOff (
            unsigned short On,
            int TimeSlotNr )
11.123.2.45 SetHeadstageSamplingActive() void SetHeadstageSamplingActive (
            bool Active,
            int TimeSlotNr )
11.123.2.46 SetHeadstageToSleep() void SetHeadstageToSleep (
            unsigned int Sleep16ms,
            int TimeSlotNr )
11.123.2.47 SetMultiHeadstageMode() void SetMultiHeadstageMode (
            bool Mode )
11.123.2.48 SetSelectedChannels() void SetSelectedChannels (
            array< BYTE >^{\wedge} channels,
            int TimeSlotNr )
```


11.123.3 Property Documentation

11.123.3.1 PulseGenerator CPulseGeneratorFunctionNet^ PulseGenerator [get]

11.123.3.2 Stimulator CW2100_StimulatorFunctionNet^ Stimulator [get]

11.124 CW2100_StimulatorFunctionNet Class Reference

Inheritance diagram for CW2100_StimulatorFunctionNet:

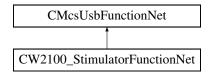

Public Member Functions

- CW2100_StimulatorFunctionNet (CMcsUsbNet[^] mcsusb)
- void SendStart (uint32_t triggermap)

Start (Trigger) the STG. The startup delay is in the range of a few ms.

void SendStop (uint32_t triggermap)

Stop some or all triggers of the STG.

- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData ^ PrepareData (int channel, array< int32_t > ^ amplitude, array< uint64_t > ^ duration, STG_DestinationEnumNet destType, uint32_t repeat)
- CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData
 PrepareDataSync (int channel, array< int32_t >^ amplitude, array< uint32_t >^ Sync, array< uint64_t >^ duration, STG_DestinationEnumNet destType, uint32_t repeat)
- void SendPreparedData (int channel, CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData[^] device_data_and_unrolled, STG_DestinationEnumNet destType)
- · void ClearChannelData (int channel)

Delete a Stimulus Pattern from STG memory

int GetDACResolution ()

Gets number of bits of the DAC resolution.

int GetTimeResolutionInNanoSeconds ()

Gets number of bits of the DAC resolution.

int GetVoltageRangeInMicroVolt (uint32_t channel)

Gets the Voltage Range of the specified channel in Microvolts.

• int GetVoltageResolutionInMicroVolt (uint32_t channel)

Gets the Voltage Resolution of the specified channel in Microvolts.

int GetCurrentRangeInNanoAmp (uint32_t channel)

Gets the Current Range of the specified channel in Nanoamps.

int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

- uint32_t GetNumberOfAnalogChannels ()
- uint32_t GetNumberOfSyncoutChannels ()
- uint32 t GetNumberOfTriggerInputs ()
- void SelectTimeSlot (int TimeSlotNr)
- int GetTimeSlot ()
- uint32_t GetStimulationPatternMemory ()

- uint32_t GetBoostPreTime ()
- uint32_t GetBoostAlwaysOnMode ()
- void SetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger number, W2100DigitalSourceEnumNet digstream source, int bitnumber offset)
- void GetDigitalStimulatorTrigger (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, [System::Runtime::InteropServices::Out]W2100DigitalSourceEnumNet% digstream_source, [System::Runtime::InteropServices::Out]int% bitnumber offset)
- void SetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number, DigitalStimulatorTriggerSlopeEnumNet slope)
- DigitalStimulatorTriggerSlopeEnumNet GetDigitalStimulatorTriggerSlope (int TimeSlotNr, DigitalStimulatorTriggerEventEnumNet trigger_event, int trigger_number)
- void StartPoll ()
- void StopPoll ()

Static Public Attributes

- static const uint32_t BOOST_BIT = (1 << 0)
- static const uint32 t GND SWITCH BIT = (1 << 1)
- static const uint32_t SYNC_BIT0 = (1 << 2)
- static const uint32_t SYNC_BIT1 = (1 << 3)

Events

OnStgPollStatus[^] PollStatusEvent

Additional Inherited Members

11.124.1 Constructor & Destructor Documentation

```
11.124.1.1 CW2100_StimulatorFunctionNet() CW2100_StimulatorFunctionNet (
CMcsUsbNet^ mcsusb )
```

11.124.2 Member Function Documentation

Delete a Stimulus Pattern from STG memory

Parameters

channel specifies the channel to clear.

```
11.124.2.2 GetBoostAlwaysOnMode() uint32_t GetBoostAlwaysOnMode ( )
```

```
11.124.2.3 GetBoostPreTime() uint32_t GetBoostPreTime ( )
```

11.124.2.4 GetCurrentRangeInNanoAmp() int GetCurrentRangeInNanoAmp (uint32_t *channel*)

Gets the Current Range of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Range of the specified channel in Nanoamps.

11.124.2.5 GetCurrentResolutionInNanoAmp() int GetCurrentResolutionInNanoAmp (uint32_t channel)

Gets the Current Resolution of the specified channel in Nanoamps.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Current Resolution of the specified channel in Nanoamps.

11.124.2.6 GetDACResolution() int GetDACResolution ()

Gets number of bits of the DAC resolution.

Returns

The DAC resolution in bits.

```
11.124.2.7 GetDigitalStimulatorTrigger() void GetDigitalStimulatorTrigger (
              int TimeSlotNr,
              {\tt DigitalStimulatorTriggerEventEnumNet}\ trigger\_event,
              int trigger_number,
              [System::Runtime::InteropServices::Out] \begin{tabular}{ll} W2100DigitalSourceEnumNet \$ & digstream\_ \hookleftarrow \end{tabular}
source,
              [System::Runtime::InteropServices::Out] int% bitnumber_offset )
\textbf{11.124.2.8} \quad \textbf{GetDigitalStimulatorTriggerSlope()} \quad \texttt{DigitalStimulatorTriggerSlopeEnumNet} \quad \texttt{GetDigital} \leftarrow \\
StimulatorTriggerSlope (
              int TimeSlotNr,
              DigitalStimulatorTriggerEventEnumNet trigger_event,
              int trigger_number )
11.124.2.9 GetNumberOfAnalogChannels() uint32_t GetNumberOfAnalogChannels ( )
11.124.2.10 GetNumberOfSyncoutChannels() uint32_t GetNumberOfSyncoutChannels ()
11.124.2.11 GetNumberOfTriggerInputs() uint32_t GetNumberOfTriggerInputs ( )
11.124.2.12 GetStimulationPatternMemory() uint32_t GetStimulationPatternMemory ( )
11.124.2.13 GetTimeResolutionInNanoSeconds() int GetTimeResolutionInNanoSeconds ( )
Gets number of bits of the DAC resolution.
Returns
     The time resolution in ns.
11.124.2.14 GetTimeSlot() int GetTimeSlot ( )
11.124.2.15 GetVoltageRangeInMicroVolt() int GetVoltageRangeInMicroVolt (
              uint32_t channel )
```

Generated by Doxygen

Gets the Voltage Range of the specified channel in Microvolts.

Parameters

channel	Channel which is queried.
---------	---------------------------

Returns

The Voltage Range of the specified channel in Microvolts.

```
11.124.2.16 GetVoltageResolutionInMicroVolt() int GetVoltageResolutionInMicroVolt ( uint32_t channel)
```

Gets the Voltage Resolution of the specified channel in Microvolts.

Parameters

channel Channel which is queried.

Returns

The Voltage Resolution of the specified channel in Microvolts.

```
11.124.2.19 SelectTimeSlot() void SelectTimeSlot ( int TimeSlotNr )
```

```
11.124.2.20 SendPreparedData() void SendPreparedData (
```

```
int channel, {\tt CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData^$ device\_data\_and\_unrolled, \\ {\tt STG\_DestinationEnumNet} \ destType \ )
```

```
11.124.2.21 SendStart() void SendStart (
uint32_t triggermap)
```

Start (Trigger) the STG. The startup delay is in the range of a few ms.

Parameters

	triggermap	A bitmap of triggers which will be started.
--	------------	---

```
11.124.2.22 SendStop() void SendStop ( uint32_t triggermap)
```

Stop some or all triggers of the STG.

Parameters

triggermap A bitmap of triggers which will be stopped.

```
11.124.2.23 SetDigitalStimulatorTrigger() void SetDigitalStimulatorTrigger (
```

```
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
W2100DigitalSourceEnumNet digstream_source,
int bitnumber_offset )
```

$\textbf{11.124.2.24} \quad \textbf{SetDigitalStimulatorTriggerSlope()} \quad \texttt{void SetDigitalStimulatorTriggerSlope} \quad \textbf{(}$

```
int TimeSlotNr,
DigitalStimulatorTriggerEventEnumNet trigger_event,
int trigger_number,
DigitalStimulatorTriggerSlopeEnumNet slope )
```

11.124.2.25 StartPoll() void StartPoll ()

```
11.124.2.26 StopPoll() void StopPoll ()
```

11.124.3 Member Data Documentation

```
11.124.3.1 BOOST_BIT const uint32_t BOOST_BIT = (1 << 0) [static]
```

```
11.124.3.2 GND_SWITCH_BIT const uint32_t GND_SWITCH_BIT = (1 << 1) [static]
```

```
11.124.3.3 SYNC_BIT0 const uint32_t SYNC_BIT0 = (1 << 2) [static]
```


```
11.124.3.4 SYNC_BIT1 const uint32_t SYNC_BIT1 = (1 << 3) [static]
```

11.124.4 Event Documentation

11.124.4.1 PollStatusEvent OnStgPollStatus^ PollStatusEvent

11.125 CW2100DacqGroupChannelSelectionNet Class Reference

Inheritance diagram for CW2100DacqGroupChannelSelectionNet:

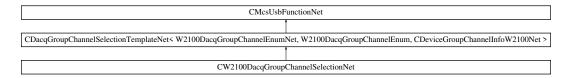

Public Member Functions

CW2100DacqGroupChannelSelectionNet (CMcsUsbNet[^] mcsusb)

Additional Inherited Members


11.125.1 Constructor & Destructor Documentation

11.125.1.1 CW2100DacqGroupChannelSelectionNet() CW2100DacqGroupChannelSelectionNet (CMcsUsbNet^ mcsusb)

11.126 CWarnerUssingDeviceNet Class Reference

CWarnerUssingDeviceNet is the class to control the Ussing device

Inheritance diagram for CWarnerUssingDeviceNet:

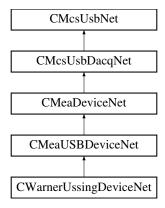

Public Member Functions

- CWarnerUssingDeviceNet ()
 Initializes a new instance of the CWarnerUssingDeviceNet class.
- virtual ∼CWarnerUssingDeviceNet ()
- !CWarnerUssingDeviceNet ()

Properties

• CWarnerUssingFunctionNet[^] WarnerUssingFunction [get]

Additional Inherited Members

11.126.1 Detailed Description

CWarnerUssingDeviceNet is the class to control the Ussing device

11.126.2 Constructor & Destructor Documentation

11.126.2.1 CWarnerUssingDeviceNet() CWarnerUssingDeviceNet ()

Initializes a new instance of the CWarnerUssingDeviceNet class.

```
11.126.2.2 ~CWarnerUssingDeviceNet() virtual ~CWarnerUssingDeviceNet () [virtual]
```

11.126.2.3 "!CWarnerUssingDeviceNet() !CWarnerUssingDeviceNet ()


11.126.3 Property Documentation

11.126.3.1 WarnerUssingFunction CWarnerUssingFunctionNet^ WarnerUssingFunction [get]

11.127 CWarnerUssingFunctionNet Class Reference

CWarnerUssingFunctionNet is the class to control the Ussing device

Inheritance diagram for CWarnerUssingFunctionNet:

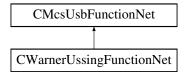

Public Member Functions

• CWarnerUssingFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] pWarner← UssingFunctionPointerContainer)

Initializes a new instance of the CWarnerUssingFunctionNet class.

- CWarnerUssingFunctionNet (CMcsUsbNet[∧] mcsusb)
- virtual ~CWarnerUssingFunctionNet ()
- !CWarnerUssingFunctionNet ()
- int32_t GetChannelsCountOfChamber (int32_t ChamberId)

gets number of channels in datastream from chamber amp with given index

int32_t GetNumberOfHardwareSlotsForChambers ()

gets number of physical hardware slots for chambers amps

int32 t GetNumberOfAvailableChambers ()

gets number of actually connected chamber amps

bool IsChamberAvailable (int32_t ChamberId)

checks whether chamber amp is connected to slot

void SetPulse (int32_t ChamberId, UssingClampModeEnumNet StgMode, int32_t NumberOfRepetitions, array< int >^ Amplitudes, array< int >^ Durations, array< int >^ PulseMarker)

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

void SetVoltageClampControllerParam_P (int32_t ChamberId, uint32_t P)

sets P value of PID controller;

void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)

sets I value of PID controller;

```
    void SetVoltageClampControllerParam_D (int32_t ChamberId, uint32_t D)

     sets D value of PID controller;

    uint32 t GetVoltageClampControllerParam P (int32 t ChamberId)

     gets P value of PID controller;

    uint32_t GetVoltageClampControllerParam_I (int32_t ChamberId)

     gets I value of PID controller;

    uint32_t GetVoltageClampControllerParam_D (int32_t ChamberId)

     gets D value of PID controller;

    void SetClampMode (int32_t Chamberld, UssingClampModeEnumNet NewClampMode)

     sets clamp mode (voltage, current or open clamp)

    UssingClampModeEnumNet GetClampMode (int32 t ChamberId)

     gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)
• bool IsInternalCalibrationFinished (int32_t ChamberId)
      when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

    int32 t GetU1Offset (int32 t ChamberId)

    int32 t GetU2Offset (int32 t ChamberId)

    int32_t GetU1Reference (int32_t Chamberld)

    int32_t GetU2Reference (int32_t Chamberld)

    int32 t GetDacZero (int32 t ChamberId)

    void SetHighCurrentMode (int32 t ChamberId)

     switch to high-current mode

    void SetLowCurrentMode (int32 t ChamberId)

     switch to low-current mode

    bool IsHighCurrentMode (int32 t ChamberId)

    uint32_t GetLowCurrentRange (int32_t Chamberld)

     query the range of the low current mode

    uint32_t GetHighCurrentRange (int32_t ChamberId)

     query the range of the high current mode

    uint32_t GetDacPampsPerDigitLowCurrentRange (int32_t ChamberId)

     get the resolution of the low current mode

    uint32 t GetDacPampsPerDigitHighCurrentRange (int32 t ChamberId)

     get the resolution of the high current mode

    uint32_t GetUnitsPerDigit (int32_t Chamberld, int32_t Channelld)

     gets amps/volts per digit for specified chamber and channel

    int32 t GetUnitExponent (int32 t Chamberld, int32 t Channelld)

     gets the unit exponent for specified chamber and channel

    UssingUnitEnumNet GetUnitName (int32_t ChamberId, int32_t ChannelId)

     gets the channel's unit name

    String \(^{\text{GetUnitDescription}}\) (int32 t Chamberld, int32 t Channelld)

     gets the description for the unit

    array< int > ^ GetAvailableChambers ()

     returns array with (zero-based) Chamberlds of all available chambers

    int32_t GetUptimeSeconds (int32_t ChamberId)

     gets the uptime in seconds

    void SetIdleModeOffset (int32_t ChamberId, UssingClampModeEnumNet ClampMode, int32_t NewIdle

  Offset)
     sets the offset (voltage or current) that will be applied when clamping is DISABLED

    int32_t GetIdleModeOffset (int32_t ChamberId, UssingClampModeEnumNet ClampMode)

     gets the offset (voltage or current) that will be applied when clamping is DISABLED

    void SetEnablePulse (int32 t Chamberld, UssingClampModeEnumNet ClampMode, bool Enable)

     enable pulse of given chamber and mode (voltage/current clamp) of this chamber
```

- bool IsPulseEnabled (int32_t Chamberld, UssingClampModeEnumNet ClampMode)
 - returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED
- void SetLiquidResistance (int32_t Chamberld, int32_t NewLiquidResistance_Ohm)

sets the resistance of the liquid

• int32 t GetLiquidResistance (int32 t Chamberld)

gets the resistance of the liquid

- int32 t GetComplianceVoltageThreshold (int32 t ChamberId)
 - returns compliance voltage threshold in uV; when Uc is above, current source is overloaded
- bool CompensateElectrodeOffset (int32_t ChamberId)

blocking call to compensate electrode offset of one chamber; returns true when successful

bool WaitForChamber (int32_t ChamberId)

blocking call that waits for chamber boot-up calibration to complete

• bool WaitForAllChambers ()

blocking call that waits for ALL chambers' boot-up calibration to complete

Additional Inherited Members

11.127.1 Detailed Description

CWarnerUssingFunctionNet is the class to control the Ussing device

11.127.2 Constructor & Destructor Documentation

```
11.127.2.1 CWarnerUssingFunctionNet() [1/2] CWarnerUssingFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ pWarnerUssingFunctionPointerContainer)
```

Initializes a new instance of the CWarnerUssingFunctionNet class.

```
11.127.2.2 CWarnerUssingFunctionNet() [2/2] CWarnerUssingFunctionNet (
CMcsUsbNet^ mcsusb )
```

```
11.127.2.3 ~CWarnerUssingFunctionNet() virtual ~CWarnerUssingFunctionNet () [virtual]
```

```
11.127.2.4 "!CWarnerUssingFunctionNet() !CWarnerUssingFunctionNet ()
```

11.127.3 Member Function Documentation

```
11.127.3.1 CompensateElectrodeOffset() bool CompensateElectrodeOffset ( int32_t ChamberId )
```

blocking call to compensate electrode offset of one chamber; returns true when successful

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

true if compensation succeeded

11.127.3.2 GetAvailableChambers() array<int> ^ GetAvailableChambers ()

returns array with (zero-based) Chamberlds of all available chambers

gets number of channels in datastream from chamber amp with given index

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

return value of zero means that amp is not placed

11.127.3.4 GetClampMode() UssingClampModeEnumNet GetClampMode (int32_t ChamberId)

gets clamp mode (voltage, current or open clamp; do not use when device is in internal calibration mode)

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

the current clamp mode

11.127.3.5 GetComplianceVoltageThreshold() int32_t GetComplianceVoltageThreshold (int32_t ChamberId)

returns compliance voltage threshold in uV; when Uc is above, current source is overloaded

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

the compliance voltage threshold in uV

get the resolution of the high current mode

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

unit: pA/digit in high current mode

11.127.3.7 GetDacPampsPerDigitLowCurrentRange() uint32_t GetDacPampsPerDigitLowCurrentRange (int32_t ChamberId)

get the resolution of the low current mode

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

pA/digit in low current mode

11.127.3.8 GetDacZero() int32_t GetDacZero (int32_t ChamberId)

 diagnostic function only -; gets real zero value of DAC in digits (0 -> neg. current; 32767 -> near zero; 65535 -> pos. current)

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

the zero value of the DAC

query the range of the high current mode

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

low current range in nA

gets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

Chamberld	index of hardware chamber slot (zero-based)
ClampMode	voltage or current clamp stimulation

Returns

unit: nA or uV

$\textbf{11.127.3.11} \quad \textbf{GetLiquidResistance()} \quad \texttt{int32_t GetLiquidResistance ()}$

int32_t ChamberId)

gets the resistance of the liquid

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

the liquid restistance in ohm

$\textbf{11.127.3.12} \quad \textbf{GetLowCurrentRange()} \quad \texttt{uint32_t} \quad \texttt{GetLowCurrentRange} \quad \textbf{(}$

int32_t ChamberId)

query the range of the low current mode

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

low current range in nA

11.127.3.13 GetNumberOfAvailableChambers() int32_t GetNumberOfAvailableChambers ()

gets number of actually connected chamber amps

Returns

the number of actually connected chambers

$\textbf{11.127.3.14} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \texttt{int32_t} \quad \texttt{GetNumberOfHardwareSlotsFor} \leftarrow \textbf{11.127.3.14} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{Int32_t} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{Int32_t} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{Int32_t} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{Int32_t} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{Int34_t} \quad \textbf{GetNumberOfHardwareSlotsForChambers()} \quad \textbf{GetNumberSlotsForChambers()} \quad \textbf{GetNumbers()}$

Chambers ()

gets number of physical hardware slots for chambers amps

Returns

the number of hardware chamber slots on the backplane

· diagnostic function only -

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

U1 offset

· diagnostic function only -

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

U1 reference

· diagnostic function only -

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

U2 offset

· diagnostic function only -

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

U2 reference

gets the description for the unit

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
Id	
Channelld	index of channel (zero-based)

Returns

the description of the unix

gets the unit exponent for specified chamber and channel

Parameters

	Chamber <i>⊷</i> Id	index of hardware chamber slot (zero-based)
Ì	Channelld	index of channel (zero-based)

Returns

example: return value -9 means that amps/volts per digit is in nano

int32_t ChannelId)

gets the channel's unit name

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	
Channelld	index of channel (zero-based)

Returns

the name of the unit

gets amps/volts per digit for specified chamber and channel

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
Id	
Channelld	index of channel (zero-based)

Returns

amps/volts per digit

```
11.127.3.23 GetUptimeSeconds() int32_t GetUptimeSeconds ( int32_t ChamberId )
```

gets the uptime in seconds

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

seconds since power-on

11.127.3.24 GetVoltageClampControllerParam_D() uint32_t GetVoltageClampControllerParam_D (int32_t ChamberId)

gets D value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)]
ld		

Returns

the D value

11.127.3.25 GetVoltageClampControllerParam_I() uint32_t GetVoltageClampControllerParam_I (int32_t ChamberId)

gets I value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

the I value

11.127.3.26 GetVoltageClampControllerParam_P() uint32_t GetVoltageClampControllerParam_P (int32_t ChamberId)

gets P value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

the P value

checks whether chamber amp is connected to slot

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	

Returns

true if the chamber is available

11.127.3.28 **IsHighCurrentMode()** bool IsHighCurrentMode (int32_t *ChamberId*)

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

true if in hight current mode

11.127.3.29 IsInternalCalibrationFinished() bool IsInternalCalibrationFinished (int32_t ChamberId)

when internal calibration is finished, values for U1,2_offset and U1,2_reference and DAC_offset are available

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

Returns

true if finished

returns true when pulse of given chamber and current mode (voltage/current clamp) of this chamber is ENABLED

Parameters

Chamberld	index of hardware chamber slot (zero-based)
ClampMode	voltage or current clamp stimulation

Returns

when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

sets clamp mode (voltage, current or open clamp)

Parameters

Chamberld	index of hardware chamber slot (zero-based)
NewClampMode	the clamp mode to use

enable pulse of given chamber and mode (voltage/current clamp) of this chamber

Parameters

Chamberld	index of hardware chamber slot (zero-based)
ClampMode	voltage or current clamp stimulation
Enable	when ENABLED, previously defined pulse pattern will be applied, otherwise the chamber current/voltage will be kept at specified offset level

```
11.127.3.33 SetHighCurrentMode() void SetHighCurrentMode ( int32_t ChamberId )
```

switch to high-current mode

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

sets the offset (voltage or current) that will be applied when clamping is DISABLED

Parameters

Chamberld	index of hardware chamber slot (zero-based)	
ClampMode	voltage or current clamp stimulation	
NewIdleOffset	unit: nA or uV	

sets the resistance of the liquid

Parameters

Chamberld	index of hardware chamber slot (zero-based)
NewLiquidResistance_Ohm	the liquid resistiance in ohm

```
11.127.3.36 SetLowCurrentMode() void SetLowCurrentMode (
    int32_t ChamberId )
```

switch to low-current mode

Chamber⊷	index of hardware chamber slot (zero-based)
ld	

```
11.127.3.37 SetPulse() void SetPulse ( int32_t ChamberId,
```

```
UssingClampModeEnumNet StgMode,
int32_t NumberOfRepetitions,
array< int >^ Amplitudes,
array< int >^ Durations,
array< int >^ PulseMarker )
```

defines stimulation pulse pattern for voltage or current stimulation; CAUTION: zero-length amplitude will be briefly applied -> choose matching to neighbour to avoid spikes

Parameters

Chamberld	index of hardware chamber slot (zero-based); send pattern to connected amp	
StgMode	voltage or current clamp stimulation	
NumberOfRepetitions	number of repetitions for pulse pattern (-1 for infinite; n means pattern is applied n+1 times)	
Amplitudes	amplitude; unit in voltage clamp: uV; unit in current clamp: nA	
Durations	duration in 100us; CAUTION: first element is applied only one; auto-loop back to second element after last one	
PulseMarker	defines values on digital channel for each step (positive: digital channel "01", neg: "10", zero: "00")	

11.127.3.38 SetVoltageClampControllerParam_D() void SetVoltageClampControllerParam_D (int32_t ChamberId, uint32_t D)

sets D value of PID controller;

Parameters

Chamber← Id	index of hardware chamber slot (zero-based)
D	useful range: 0700

11.127.3.39 SetVoltageClampControllerParam_I() void SetVoltageClampControllerParam_I (int32_t ChamberId, uint32_t I)

sets I value of PID controller;

	index of hardware chamber slot (zero-based)
ld	
1	useful range: 80000120000

11.127.3.40 SetVoltageClampControllerParam_P() void SetVoltageClampControllerParam_P (

```
int32_t ChamberId,
uint32_t P )
```

sets P value of PID controller;

Parameters

Chamber⊷	index of hardware chamber slot (zero-based)
ld	
Р	useful value: 130000

11.127.3.41 WaitForAllChambers() bool WaitForAllChambers ()

blocking call that waits for ALL chambers' boot-up calibration to complete

Returns

returns false when at least one chamber's calibration fails (e.g. timeout...)

blocking call that waits for chamber boot-up calibration to complete

Parameters

Chamber←	index of hardware chamber slot (zero-based)
ld	


Returns

returns false when calibration fails (e.g. timeout...)

11.128 CWarnerValveControllerDeviceNet Class Reference

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

Inheritance diagram for CWarnerValveControllerDeviceNet:

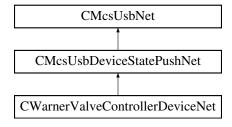

Public Member Functions

- delegate void OnGetValveActive (uint16 t valve, int valveActive)
- delegate void OnGetValveManualState (uint16_t valve, int32_t valveManualState)
- delegate void OnGetValveManualGroup (uint16_t valve, int32_t valveManualGroup)
- delegate void OnGetValveMode (uint16 t valve, WvcValveModeEnumNet ValveMode)
- delegate void OnGetAnalogThresholdLow (uint16_t valve, int32_t threshold)
- delegate void OnGetAnalogThresholdHigh (uint16_t valve, int32_t threshold)
- delegate void OnGetDigitalPortDirection (uint16_t port, PortDirectionEnumNet direction)
- delegate void OnlsValveDigitalInInverted (uint16_t valve, bool isInverted)
- delegate void OnGetValveDigitalInPort (uint16_t valve, uint32_t digitalInPort)
- delegate void OnlsDigitalOutPortInverted (uint16_t digitalOutPort, bool isInverted)
- delegate void OnGetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)
- delegate void OnlsValveOpen (uint16_t valve, bool valveOpen)
- delegate void OnlsValveOpenInDigitalMode (uint16 t valve, bool valveOpen)
- delegate void OnlsValveOpenInAnalogMode (uint16_t valve, bool valveOpen)
- delegate void OnGetAnalogVoltage (int32 t voltage)
- delegate void OnTableEntryChanged (uint16 t tableNumber)
- delegate void OnGetTableNamebyIndex (uint16 t tableNumber, String[^] tableName)
- delegate void OnGetActiveRunningTableNumber (uint32_t tableNumber)
- delegate void OnGetCurrentNumberOfValves (int32_t numberOfValves)
- delegate void OnGetValveBoardRevision (uint32_t revision)
- delegate void OnGetValveLedOn (bool ledon)
- delegate void OnGetDisplayMode (WvcDisplayModeEnumNet DisplayMode)
- CWarnerValveControllerDeviceNet ()

Initializes a new instance of the CWarnerValveControllerDeviceNet class.

- virtual ~CWarnerValveControllerDeviceNet ()
- !CWarnerValveControllerDeviceNet ()
- int GetValveActive (uint16_t valve)

Gets the valve active/inactive state

void SetValveActive (uint16_t valve, int valveActive)

Sets the valve active/inactive state

uint32_t GetValvesActiveMap ()

Gets the valves active/inactive states

void SetValvesActiveMap (uint32_t valvesActive)

Sets the valve active/inactive state

int32_t GetValveManualState (uint16_t valve)

Gets the valve manual on/off state

void SetValveManualState (uint16_t valve, int32_t valveManualState)

Sets the valve manual on/off state

uint32 t GetValvesManualStateMap ()

Gets the valves manual on/off states

void SetValvesManualStateMap (uint32 t valveaManualState)

Sets the valve manual on/off state

int32_t GetValveManualGroup (uint16_t valve)

Gets the valve manual group

void SetValveManualGroup (uint16 t valve, int32 t valveManualGroup)

Sets the valve manual group

WvcValveModeEnumNet GetValveMode (uint16_t valve)

Reads the valve mode

void SetValveMode (uint16 t valve, WvcValveModeEnumNet ValveMode)

Writes the valve mode

• int32_t GetAnalogThresholdLow (uint16_t valve)

Gets the lower threshold for the analog in port per valve

void SetAnalogThresholdLow (uint16_t valve, int32_t threshold)

Sets the lower threshold for the analog in port per valve

int32 t GetAnalogThresholdHigh (uint16 t valve)

Gets the upper threshold for the analog in port per valve

void SetAnalogThresholdHigh (uint16_t valve, int32_t threshold)

Sets the upper threshold for the analog in port per valve

PortDirectionEnumNet GetDigitalPortDirection (uint16 t port)

Gets the direction of a digital port

void SetDigitalPortDirection (uint16_t port, PortDirectionEnumNet direction)

Sets the direction of a digital port

bool IsValveDigitalInInverted (uint16_t valve)

Is digital in inverted

void SetValveDigitalInInvert (uint16_t valve, bool isInverted)

Invert digital in

uint32_t GetValveDigitalInPort (uint16_t valve)

Gets the number of the digital in port which is mapped to a valve

• void SetValveDigitalInPort (uint16_t valve, uint32_t digitalInPort)

Map a digital in port to a valve

bool IsDigitalOutPortInverted (uint16_t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

void SetDigitalOutPortInvert (uint16_t digitalOutPort, bool isInverted)

Map a valve to a digital out port

uint32_t GetDigitalOutPortValve (uint16_t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

void SetDigitalOutPortValve (uint16_t digitalOutPort, uint32_t valve)

Map a valve to a digital out port

void SetDefault ()

Sets the settings of the valve controller to default

bool IsValveOpen (uint16_t valve)

Is valve open

bool IsValveOpenInDigitalMode (uint16_t valve)

True, if the valve would be open when the device is in digital mode

bool IsValveOpenInAnalogMode (uint16_t valve)

True, if the valve would be open when the device is in analog mode

• int32_t GetAnalogVoltage ()

Reads the voltage on the analog in port

void GetValveTableEntry (uint16_t valve, uint16_t index, [System::Runtime::InteropServices::Out]uint32_t% duration, [System::Runtime::InteropServices::Out]bool% state)

Read an entry from the valve protocol table

• void SetValveTableEntry (uint16 t valve, uint16 t index, uint32 t duration, bool state)

Write an entry to the valve protocol table

void ClearValveTable (uint16_t valve)

Clear the valve protocol table

void LoadValveTable ()

Load the current table from permanent memory

void StoreValveTable ()

Store the current table in permanent memory

String \(^\) GetTableNamebyIndex (uint16 t tableNumber)

Get the name of a protocol table

• String ^ GetTableName ()

Get the name of the current protocol table

void SetTableName (String[^] tableName)

Set the name of the current protocol table

• uint32 t GetActiveRunningTableNumber ()

Gets the number of the table that is active for running

void SetActiveRunningTableNumber (uint32_t tableNumber)

Sets the number of the tanle that is active for running

uint32 t GetCurrentEditTableNumber ()

Gets the number of the table that is current for editing

void SetCurrentEditTableNumber (uint32_t tableNumber)

Sets the number of the table that is current for editing

void ClearTableName ()

Clear the name of current protocol table

• void SetTableStep (uint16_t valve, int32_t steps)

Skips the table protocol for a valve by steps

void SetTableStepAll (int32 t steps)

Skips the table protocol for all valves by steps

int32 t GetTotalNumberOfValves ()

Get the total number of valves in the system

int32_t GetTotalNumberOfDigitalPorts ()

Get the total number of digital ports in the system

int32_t GetTotalTableSize ()

Get the total table size in the system

int32_t GetTotalNumberOfTables ()

Get the total number of tables in the system

int32_t GetCurrentNumberOfValves ()

Get the current number of valves connected to the system

uint32_t GetValveBoardRevision ()

Gets the revision code of the valve board

• bool GetValveLedOn ()

Gets the LED state of the valve board

void SetValveLedOn (bool ledon)

Gets the LED state of the valve board

WvcDisplayModeEnumNet GetDisplayMode ()

Reads the display mode

void SetDisplayMode (WvcDisplayModeEnumNet DisplayMode, int32_t lockTimeMs)

Writes the display mode

String \(^\) GetValveBoardRevisionString ()

Gets the revision name of the valve board

Gets the valve currents

void SetValveCurrent (int16_t switch_current, int16_t hold_current)

Sets the valve currents different from the default

Events

```
• OnGetValveActive^ GetValveActiveEvent [add, remove, raise]
```

Event fires when the valve state for the valve number has changed

• OnGetValveManualState^ GetValveManualStateEvent [add, remove, raise]

Event fires when the manual valve state for the valve number has changed

OnGetValveManualGroup[^] GetValveManualGroupEvent [add, remove, raise]

Event fires when the manual valve group for the valve number has changed

• OnGetValveMode GetValveModeEvent [add, remove, raise]

Event fires when the valve mode for the valve number has changed

• OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent [add, remove, raise]

Event fires when the threshold in mV for the valve number has changed

OnGetAnalogThresholdHigh[^] GetAnalogThresholdHighEvent [add, remove, raise]

Event fires when the threshold in mV for the valve number has changed

• OnGetDigitalPortDirection GetDigitalPortDirectionEvent [add, remove, raise]

Event fires when the direction for the port number has changed

• OnlsValveDigitalInInverted^ IsValveDigitalInInvertedEvent [add, remove, raise]

Event fires when is inverted for the valve number has changed

• OnGetValveDigitalInPort^ GetValveDigitalInPortEvent [add, remove, raise]

Event fires when the digital in port for the valve number has changed

• OnlsDigitalOutPortInverted^ IsDigitalOutPortInvertedEvent [add, remove, raise]

Event fires when is inverted for the digital out port has changed

• OnGetDigitalOutPortValve GetDigitalOutPortValveEvent [add, remove, raise]

Event fires when the valve number for the digital out port has changed

• OnlsValveOpen^ IsValveOpenEvent [add, remove, raise]

Event fires when is open for the valve number has changed

OnlsValveOpenInDigitalMode^ IsValveOpenInDigitalModeEvent [add, remove, raise]

Event fires when is open for the valve number has changed

• OnlsValveOpenInAnalogMode^ IsValveOpenInAnalogModeEvent [add, remove, raise]

Event fires when is open for the valve number has changed

• OnGetAnalogVoltage GetAnalogVoltageEvent [add, remove, raise]

Event fires when the voltage in mV has changed

OnTableEntryChanged^ TableEntryChangedEvent [add, remove, raise]

Event fires when an entry of a table changed

• OnGetTableNamebyIndex^ GetTableNamebyIndexEvent [add, remove, raise]

Event fires when the name of the table for the table number has changed

OnGetActiveRunningTableNumber[^] GetActiveRunningTableNumberEvent [add, remove, raise]

Event fires when the table number has changed

• OnGetCurrentNumberOfValves^ GetCurrentNumberOfValvesEvent [add, remove, raise]

Event fires when the number of valves has changed

• OnGetValveBoardRevision GetValveBoardRevisionEvent [add, remove, raise]

Event fires when the revision code has changed

• OnGetValveLedOn^ GetValveLedOnEvent [add, remove, raise]

Event fires when the LED state has changed

• OnGetDisplayMode^ GetDisplayModeEvent [add, remove, raise]

Event fires when the display mode has changed

Additional Inherited Members

11.128.1 Detailed Description

CWarnerValveControllerDeviceNet is the class to access the Warner Valve Controller

11.128.2 Constructor & Destructor Documentation

```
11.128.2.1 CWarnerValveControllerDeviceNet() CWarnerValveControllerDeviceNet ()
```

Initializes a new instance of the CWarnerValveControllerDeviceNet class.

```
11.128.2.2 ~CWarnerValveControllerDeviceNet() virtual ~CWarnerValveControllerDeviceNet () [virtual]
```

```
11.128.2.3 "!CWarnerValveControllerDeviceNet() !CWarnerValveControllerDeviceNet ()
```

11.128.3 Member Function Documentation

```
11.128.3.1 ClearTableName() void ClearTableName ( )
```

Clear the name of current protocol table

```
11.128.3.2 ClearValveTable() void ClearValveTable ( uint16_t valve )
```

Clear the valve protocol table

Parameters

valve The valve number

11.128.3.3 GetActiveRunningTableNumber() uint32_t GetActiveRunningTableNumber ()

Gets the number of the table that is active for running

Returns

The table number

11.128.3.4	GetAnalogThresholdHigh()	<pre>int32_t GetAnalogThresholdHigh</pre>	(
	uint16_t <i>valve</i>)		

Gets the upper threshold for the analog in port per valve

Parameters

Returns

The threshold in mV

11.128.3.5 GetAnalogThresholdLow() int32_t GetAnalogThresholdLow (uint16_t valve)

Gets the lower threshold for the analog in port per valve

Parameters

valve The valve number

Returns

The threshold in mV

11.128.3.6 GetAnalogVoltage() int32_t GetAnalogVoltage ()

Reads the voltage on the analog in port

Returns

The voltage in mV

11.128.3.7 GetCurrentEditTableNumber() uint32_t GetCurrentEditTableNumber ()

Gets the number of the table that is current for editing

Returns

The table number

11.128.3.8 GetCurrentNumberOfValves() int32_t GetCurrentNumberOfValves ()

Get the current number of valves connected to the system

Returns

The number of valves

```
11.128.3.9 GetDigitalOutPortValve() uint32_t GetDigitalOutPortValve ( uint16_t digitalOutPort )
```

Gets the number of the valve which is mapped to a digital out port

Parameters

digitalOutPort	The digital out port
----------------	----------------------

Returns

The valve number

11.128.3.10 GetDigitalPortDirection() PortDirectionEnumNet GetDigitalPortDirection (uint16_t port)

Gets the direction of a digital port

Parameters

```
port The port number
```

Returns

the direction

11.128.3.11 GetDisplayMode() WvcDisplayModeEnumNet GetDisplayMode ()

Reads the display mode

Returns

The display mode

11.128.3.12 GetTableName()	String ^	GetTableName	()
----------------------------	----------	--------------	---	---

Get the name of the current protocol table

Returns

The name of the table

```
11.128.3.13 GetTableNamebyIndex() String ^{\land} GetTableNamebyIndex ( uint16_t tableNumber )
```

Get the name of a protocol table

Parameters

tableNumber The table number	
------------------------------	--

Returns

The name of the table

11.128.3.14 GetTotalNumberOfDigitalPorts() int32_t GetTotalNumberOfDigitalPorts ()

Get the total number of digital ports in the system

Returns

The number of digital ports

11.128.3.15 GetTotalNumberOfTables() int32_t GetTotalNumberOfTables ()

Get the total number of tables in the system

Returns

The number of tables

11.128.3.16 GetTotalNumberOfValves() int32_t GetTotalNumberOfValves ()

Get the total number of valves in the system

Returns

The number of valves

11.128.3.17 GetTotalTableSize() int32_t GetTotalTableSize ()

Get the total table size in the system

Returns

The table size

11.128.3.18 GetValveActive() int GetValveActive (uint16_t valve)

Gets the valve active/inactive state

Parameters

valve	The valve number
-------	------------------

Returns

The valve state

11.128.3.19 GetValveBoardRevision() uint32_t GetValveBoardRevision ()

Gets the revision code of the valve board

Returns

The revision code

11.128.3.20 GetValveBoardRevisionString() String ^ GetValveBoardRevisionString ()

Gets the revision name of the valve board

Returns

The revision name

Gets the valve currents

Parameters

switch_current	The switch current (in DAC units)
hold_current	The hold current (in DAC units)

```
11.128.3.22 GetValveDigitalInPort() uint32_t GetValveDigitalInPort ( uint16_t valve )
```

Gets the number of the digital in port which is mapped to a valve

Parameters

Returns

The digital in port

$\textbf{11.128.3.23} \quad \textbf{GetValveLedOn()} \quad \texttt{bool GetValveLedOn ()} \\$

Gets the LED state of the valve board

Returns

The LED state

```
11.128.3.24 GetValveManualGroup() int32_t GetValveManualGroup ( uint16_t valve )
```

Gets the valve manual group

valve The valve number

Returns

The manual valve group

```
11.128.3.25 GetValveManualState() int32_t GetValveManualState ( uint16_t valve )
```

Gets the valve manual on/off state

Parameters

valve	The valve number
-------	------------------

Returns

The manual valve state

```
11.128.3.26 GetValveMode() WvcValveModeEnumNet GetValveMode ( uint16_t valve )
```

Reads the valve mode

Parameters

valve The valve number

Returns

The valve mode

11.128.3.27 GetValvesActiveMap() uint32_t GetValvesActiveMap ()

Gets the valves active/inactive states

Returns

The valves states

11.128.3.28 GetValvesManualStateMap() uint32_t GetValvesManualStateMap ()

Gets the valves manual on/off states

Returns

The manual valves states

```
\textbf{11.128.3.29} \quad \textbf{GetValveTableEntry()} \quad \texttt{void GetValveTableEntry} \ \ (
```

```
uint16_t valve,
uint16_t index,
[System::Runtime::InteropServices::Out] uint32_t% duration,
[System::Runtime::InteropServices::Out] bool% state )
```

Read an entry from the valve protocol table

Parameters

valve	The valve number	
index	The index in the table	
duration	the duration in ms	
state	the state	

11.128.3.30 IsDigitalOutPortInverted() bool IsDigitalOutPortInverted (uint16_t digitalOutPort)

Gets the number of the valve which is mapped to a digital out port

Parameters

digitalOutPort	The digital out port

Returns

is inverted

11.128.3.31 IsValveDigitalInInverted() bool IsValveDigitalInInverted (uint16_t *valve*)

Is digital in inverted

e The valve	e The valve number
e The valve	$e \mid$ The valve numbe

Returns

is inverted

```
11.128.3.32 IsValveOpen() bool IsValveOpen ( uint16_t valve )
```

Is valve open

Parameters

valve The valve nun

Returns

is open

11.128.3.33 IsValveOpenInAnalogMode() bool IsValveOpenInAnalogMode (uint16_t *valve*)

True, if the valve would be open when the device is in analog mode

Parameters

valve The valve num	ber
---------------------	-----

Returns

is open

11.128.3.34 IsValveOpenInDigitalMode() bool IsValveOpenInDigitalMode (uint16_t valve)

True, if the valve would be open when the device is in digital mode

Parameters

valve The valve number

Returns

is open

```
11.128.3.35 LoadValveTable() void LoadValveTable ()
Load the current table from permanent memory
11.128.3.36 OnGetActiveRunningTableNumber() delegate void OnGetActiveRunningTableNumber (
             uint32_t tableNumber )
11.128.3.37 OnGetAnalogThresholdHigh() delegate void OnGetAnalogThresholdHigh (
             uint16_t valve,
             int32_t threshold )
11.128.3.38 OnGetAnalogThresholdLow() delegate void OnGetAnalogThresholdLow (
             uint16_t valve,
             int32_t threshold )
11.128.3.39 OnGetAnalogVoltage() delegate void OnGetAnalogVoltage (
             int32_t voltage )
11.128.3.40 OnGetCurrentNumberOfValves() delegate void OnGetCurrentNumberOfValves (
             int32_t numberOfValves )
\textbf{11.128.3.41} \quad \textbf{OnGetDigitalOutPortValve()} \quad \texttt{delegate void OnGetDigitalOutPortValve ()}
             uint16_t digitalOutPort,
             uint32_t valve )
11.128.3.42 OnGetDigitalPortDirection() delegate void OnGetDigitalPortDirection (
             uint16_t port,
             PortDirectionEnumNet direction )
11.128.3.43 OnGetDisplayMode() delegate void OnGetDisplayMode (
             WvcDisplayModeEnumNet DisplayMode )
```

```
11.128.3.44 OnGetTableNamebyIndex() delegate void OnGetTableNamebyIndex (
             uint16_t tableNumber,
             String^{\wedge} tableName )
11.128.3.45 OnGetValveActive() delegate void OnGetValveActive (
             uint16_t valve,
             int valveActive )
11.128.3.46 OnGetValveBoardRevision() delegate void OnGetValveBoardRevision (
             uint32_t revision )
11.128.3.47 OnGetValveDigitalInPort() delegate void OnGetValveDigitalInPort (
             uint16_t valve,
             uint32_t digitalInPort )
11.128.3.48 OnGetValveLedOn() delegate void OnGetValveLedOn (
             bool ledon )
11.128.3.49 OnGetValveManualGroup() delegate void OnGetValveManualGroup (
             uint16_t valve,
             int32_t valveManualGroup )
11.128.3.50 OnGetValveManualState() delegate void OnGetValveManualState (
             uint16_t valve,
             int32_t valveManualState )
11.128.3.51 OnGetValveMode() delegate void OnGetValveMode (
             uint16_t valve,
             WvcValveModeEnumNet ValveMode )
```

```
11.128.3.52 OnlsDigitalOutPortInverted() delegate void OnlsDigitalOutPortInverted (
               uint16_t digitalOutPort,
               bool isInverted )
\textbf{11.128.3.53} \quad \textbf{OnlsValveDigitalInInverted()} \quad \texttt{delegate void OnlsValveDigitalInInverted ()}
               uint16_t valve,
               bool is Inverted )
11.128.3.54 OnlsValveOpen() delegate void OnIsValveOpen (
              uint16_t valve,
               bool valveOpen )
\textbf{11.128.3.55} \quad \textbf{OnlsValveOpenInAnalogMode()} \quad \texttt{delegate void OnlsValveOpenInAnalogMode ()}
               uint16_t valve,
               bool valveOpen )
11.128.3.56 OnlsValveOpenInDigitalMode() delegate void OnlsValveOpenInDigitalMode (
              uint16_t valve,
              bool valveOpen )
11.128.3.57 OnTableEntryChanged() delegate void OnTableEntryChanged (
               uint16_t tableNumber )
11.128.3.58 SetActiveRunningTableNumber() void SetActiveRunningTableNumber (
               uint32_t tableNumber )
Sets the number of the tanle that is active for running
Parameters
 tableNumber
                The table number
\textbf{11.128.3.59} \quad \textbf{SetAnalogThresholdHigh()} \quad \texttt{void SetAnalogThresholdHigh ()}
              uint16_t valve,
```

int32_t threshold)

Sets the upper threshold for the analog in port per valve

Parameters

valve	The valve number
threshold	The threshold in mV

Sets the lower threshold for the analog in port per valve

Parameters

valve	The valve number
threshold	The threshold in mV

11.128.3.61 SetCurrentEditTableNumber() void SetCurrentEditTableNumber (uint32_t tableNumber)

Sets the number of the table that is current for editing

Parameters

tableNumber	The table number

$\textbf{11.128.3.62} \quad \textbf{SetDefault()} \quad \texttt{void SetDefault ()}$

Sets the settings of the valve controller to default

Map a valve to a digital out port

digitalOutPort	The digital out port
isInverted	True if digital out is to be inverted

```
11.128.3.64 SetDigitalOutPortValve() void SetDigitalOutPortValve ( uint16_t digitalOutPort, uint32_t valve )
```

Map a valve to a digital out port

Parameters

digitalOutPort	The digital out port
valve	The valve number

Sets the direction of a digital port

Parameters

port	The port number
direction	the direction

Writes the display mode

Parameters

DisplayMode	The display mode
lockTimeMs	Locks the display for ms

```
11.128.3.67 SetTableName() void SetTableName ( String<sup>^</sup> tableName )
```

Set the name of the current protocol table

|--|

Skips the table protocol for a valve by steps

Parameters

valve	The valve number
steps	Number of steps

Skips the table protocol for all valves by steps

Parameters

steps	Number of steps
-------	-----------------

Sets the valve active/inactive state

Parameters

valve	The valve number
valveActive	The valve state

Sets the valve currents different from the default

switch_current	The switch current (in DAC units); -1 sets the device default current	
hold_current	The hold current (in DAC units); -1 sets the device default current	

11.128.3.72 SetValveDigitalInInvert() void SetValveDigitalInInvert (

```
uint16_t valve,
bool isInverted )
```

Invert digital in

Parameters

valve	The valve number
isInverted	True if digital in is to be inverted

$\textbf{11.128.3.73} \quad \textbf{SetValveDigitalInPort()} \quad \texttt{void SetValveDigitalInPort} \quad \textbf{(}$

```
uint16_t valve,
uint32_t digitalInPort )
```

Map a digital in port to a valve

Parameters

valve	The valve number
digitalInPort	The digital in port

11.128.3.74 SetValveLedOn() void SetValveLedOn (

bool ledon)

Gets the LED state of the valve board

Parameters

ledon The LED state

11.128.3.75 SetValveManualGroup() void SetValveManualGroup (

```
uint16_t valve,
int32_t valveManualGroup )
```

Sets the valve manual group

valve	The valve number
valveManualGroup	The manual valve group

Sets the valve manual on/off state

Parameters

valve	The valve number
valveManualState	The manual valve state

Writes the valve mode

Parameters

valve	The valve number
ValveMode	The valve mode

11.128.3.78 SetValvesActiveMap() void SetValvesActiveMap (uint32_t valvesActive)

Sets the valve active/inactive state

Parameters

valvesActive	The valves states

11.128.3.79 SetValvesManualStateMap() void SetValvesManualStateMap (uint32_t valveaManualState)

Sets the valve manual on/off state

valveaManualState	The manual valves states
-------------------	--------------------------

11.128.3.80 SetValveTableEntry() void SetValveTableEntry (

```
uint16_t valve,
uint16_t index,
uint32_t duration,
bool state )
```

Write an entry to the valve protocol table

Parameters

valve	The valve number
index	The index in the table
duration	the duration in ms
state	the state

11.128.3.81 StoreValveTable() void StoreValveTable ()

Store the current table in permanent memory

11.128.4 Event Documentation

11.128.4.1 GetActiveRunningTableNumberEvent OnGetActiveRunningTableNumber^ GetActiveRunning← TableNumberEvent [add], [remove], [raise]

Event fires when the table number has changed

11.128.4.2 GetAnalogThresholdHighEvent OnGetAnalogThresholdHigh^ GetAnalogThresholdHighEvent [add], [remove], [raise]

Event fires when the threshold in mV for the valve number has changed

11.128.4.3 GetAnalogThresholdLowEvent OnGetAnalogThresholdLow^ GetAnalogThresholdLowEvent [add], [remove], [raise]

Event fires when the threshold in mV for the valve number has changed

11.128.4.4 GetAnalogVoltageEvent OnGetAnalogVoltage^ GetAnalogVoltageEvent [add], [remove], [raise]

Event fires when the voltage in mV has changed

11.128.4.5 GetCurrentNumberOfValvesEvent OnGetCurrentNumberOfValves^ GetCurrentNumberOf↔ ValvesEvent [add], [remove], [raise]

Event fires when the number of valves has changed

11.128.4.6 GetDigitalOutPortValveEvent OnGetDigitalOutPortValve^ GetDigitalOutPortValveEvent [add], [remove], [raise]

Event fires when the valve number for the digital out port has changed

11.128.4.7 GetDigitalPortDirectionEvent OnGetDigitalPortDirection^ GetDigitalPortDirectionEvent [add], [remove], [raise]

Event fires when the direction for the port number has changed

11.128.4.8 **GetDisplayModeEvent** OnGetDisplayMode^ GetDisplayModeEvent [add], [remove], [raise]

Event fires when the display mode has changed

11.128.4.9 GetTableNamebyIndexEvent OnGetTableNamebyIndex^ GetTableNamebyIndexEvent [add], [remove], [raise]

Event fires when the name of the table for the table number has changed

11.128.4.10 GetValveActiveEvent OnGetValveActive^ GetValveActiveEvent [add], [remove], [raise]

Event fires when the valve state for the valve number has changed

11.128.4.11 GetValveBoardRevisionEvent OnGetValveBoardRevision^ GetValveBoardRevisionEvent [add], [remove], [raise]

Event fires when the revision code has changed

11.128.4.12 GetValveDigitalInPortEvent OnGetValveDigitalInPort^ GetValveDigitalInPortEvent [add], [remove], [raise]

Event fires when the digital in port for the valve number has changed

 $\textbf{11.128.4.13} \quad \textbf{GetValveLedOnEvent} \quad \texttt{OnGetValveLedOn}^{\land} \quad \texttt{GetValveLedOnEvent} \quad \texttt{[add], [remove], [raise]}$

Event fires when the LED state has changed

11.128.4.14 GetValveManualGroupEvent OnGetValveManualGroup^ GetValveManualGroupEvent [add], [remove], [raise]

Event fires when the manual valve group for the valve number has changed

11.128.4.15 GetValveManualStateEvent OnGetValveManualState^ GetValveManualStateEvent [add], [remove], [raise]

Event fires when the manual valve state for the valve number has changed

 $\textbf{11.128.4.16} \quad \textbf{GetValveModeEvent} \quad \texttt{OnGetValveMode}^{\wedge} \quad \texttt{GetValveModeEvent} \quad \texttt{[add], [remove], [raise]}$

Event fires when the valve mode for the valve number has changed

11.128.4.17 IsDigitalOutPortInvertedEvent OnIsDigitalOutPortInverted^ IsDigitalOutPortInverted← Event [add], [remove], [raise]

Event fires when is inverted for the digital out port has changed

11.128.4.18 IsValveDigitalInInvertedEvent OnIsValveDigitalInInverted[∧] IsValveDigitalInInverted← Event [add], [remove], [raise]

Event fires when is inverted for the valve number has changed

11.128.4.19 IsValveOpenEvent OnIsValveOpen^ IsValveOpenEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.128.4.20 IsValveOpenInAnalogModeEvent OnIsValveOpenInAnalogMode^ IsValveOpenInAnalogMode← Event [add], [remove], [raise]

Event fires when is open for the valve number has changed

11.128.4.21 IsValveOpenInDigitalModeEvent OnIsValveOpenInDigitalMode^ IsValveOpenInDigital← ModeEvent [add], [remove], [raise]

Event fires when is open for the valve number has changed


11.128.4.22 TableEntryChangedEvent OnTableEntryChanged^ TableEntryChangedEvent [add], [remove], [raise]

Event fires when an entry of a table changed

11.129 CWarnerValveControllerDeviceTesterFunctionNet Class Reference

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

 $Inheritance\ diagram\ for\ CWarner Valve Controller Device Tester Function Net:$

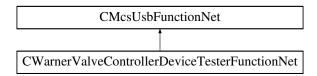

Public Member Functions

CWarnerValveControllerDeviceTesterFunctionNet (CMcsUsbNet[^] mcsusb, CMcsUsbFunctionPointerContainer[^] pWarnerValveControllerDeviceTesterFunctionPointerContainer)

Initializes a new instance of the CWarnerValveControllerDeviceTesterFunctionNet class.

- CWarnerValveControllerDeviceTesterFunctionNet (CMcsUsbNet[^] mcsusb)
- virtual ~CWarnerValveControllerDeviceTesterFunctionNet ()
- !CWarnerValveControllerDeviceTesterFunctionNet ()
- void SetADC (uint32_t onoff)

Sets the ADC port of the tester

• uint32 t GetSync ()

Gets the output from the sync port

void SetTrigger (uint32 t trigger)

Sets the input to the trigger port

void SetTriggerSyncDirection (uint32 t direction)

Sets the direction of the trigger/sync test port

· uint32 t GetIO ()

Gets the output from the io ports

void SetIO (uint32_t io)

Sets the input to the io ports

· void SetIODirection (int32 t direction)

Sets the direction of the IO test ports

Additional Inherited Members

11.129.1 Detailed Description

CWarnerValveControllerDeviceTesterFunctionNet is the class to access the functions for the Warner Valve Controller Device Tester

11.129.2 Constructor & Destructor Documentation

```
11.129.2.1 CWarnerValveControllerDeviceTesterFunctionNet() [1/2] CWarnerValveControllerDeviceTesterFunctionNet (

CMcsUsbNet^ mcsusb,

CMcsUsbFunctionPointerContainer^ pWarnerValveControllerDeviceTesterFunction↔

PointerContainer )
```

Initializes a new instance of the CWarnerValveControllerDeviceTesterFunctionNet class.

```
11.129.2.2 CWarnerValveControllerDeviceTesterFunctionNet() [2/2] CWarnerValveControllerDeviceTesterFunctionNet (

CMcsUsbNet^ mcsusb )
```

```
11.129.2.3 ~CWarnerValveControllerDeviceTesterFunctionNet() virtual ~CWarnerValveControllerDeviceTesterFunction() [virtual]

11.129.2.4 "!CWarnerValveControllerDeviceTesterFunctionNet() !CWarnerValveControllerDeviceTesterFunctionNet()

11.129.3.1 GetIO() uint32_t GetIO()

Gets the output from the io ports

Returns

The manual valves states
```

11.129.3.2 GetSync() uint32_t GetSync ()

Gets the output from the sync port

Returns

The sync state

```
11.129.3.3 SetADC() void SetADC ( uint32_t onoff )
```

Sets the ADC port of the tester

Parameters

onoff The port state

```
11.129.3.4 SetIO() void SetIO ( uint32_t io )
```

Sets the input to the io ports

Parameters

io The manual valves states

11.129.3.5 SetIODirection() void SetIODirection (int32_t direction)

Sets the direction of the IO test ports

Parameters

direction The 16bit direction map: 1=IN 0=OUT

Sets the input to the trigger port

Parameters

trigger The trigger state

11.129.3.7 SetTriggerSyncDirection() void SetTriggerSyncDirection (uint32_t direction)


Sets the direction of the trigger/sync test port

Parameters

direction The direction: 1=IN 0=OUT

11.130 CWClassicFunctionNet Class Reference

Inheritance diagram for CWClassicFunctionNet:

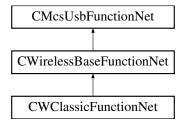

Public Member Functions

- CWClassicFunctionNet (CMcsUsbNet[∧] mcsusb, CMcsUsbFunctionPointerContainer[∧] wClassicFuntion
 — PointerContainer)
- CWClassicFunctionNet (CMcsUsbNet[^] mcsusb)
- uint32 t ResetChannelmap (unsigned int virtualDevice)
- uint32_t SetChannelmap (unsigned char position, unsigned char channel, unsigned int Device)
- void SetHWSelectedChannels (array< bool >^ channels, unsigned int Device)
- void SetRFLostBehaviour (uint8_t stoponfailure, unsigned int Device)
- · void SetHeadstageOnOff (uint16 t onoff)
- USHORT GetHeadstageOnOff ()
- void SetRFFrequencyHeadstage (uint8_t receiver_nb, unsigned short frequency)
- unsigned short GetRFFrequencyHeadstage (uint8_t receiver_nb)
- void SetRFFrequencyReceiver (uint8 t receiver nb, uint8 t configuration, unsigned short frequency)
- void SetRFFrequencyReceiverEeprom (uint8 t receiver nb, uint8 t configuration, unsigned short frequency)
- unsigned short GetRFFrequencyReceiver (uint8_t receiver_nb, uint8_t configuration)
- void SetSerialNumberHeadstage (unsigned short number)
- unsigned short GetSerialNumberHeadstage ()
- void SetSelectedHeadstage (uint8_t number)
- uint8_t GetSelectedHeadstage ()
- void ScanForHeadstages ()
- uint8_t GetScanHeadstagesResult (int max_wait_for_ms)
- void SetFilterParametersHeadstage (unsigned short index, array< int >^ buffer)
- array< int > ^ GetFilterParametersHeadstage (unsigned short index)
- bool GetHasRedLedHeadstage ()
- void SetHasChecksum (unsigned int has, unsigned int Device)
- unsigned int GetHasChecksum (unsigned int Device)
- · void SetResetFilter (unsigned int reset, unsigned int Device)
- unsigned int GetResetFilter (unsigned int Device)
- · void SetWPAType (unsigned short type, unsigned int Device)
- unsigned short GetWPAType (unsigned int Device)
- void SetWPADebugMode (unsigned int mode, unsigned int Device)
- unsigned int GetWPADebugMode (unsigned int Device)
- void SetRFPower (unsigned short power)
- unsigned short GetRFPower ()
- unsigned int GetRFConnectionStatus ()

Additional Inherited Members

11.130.1 Constructor & Destructor Documentation

```
11.130.1.1 CWClassicFunctionNet() [1/2] CWClassicFunctionNet (
              CMcsUsbNet^ mcsusb,
              {\tt CMcsUsbFunctionPointerContainer}^{\land} \ \textit{wClassicFuntionPointerContainer} \ )
11.130.1.2 CWClassicFunctionNet() [2/2] CWClassicFunctionNet (
              CMcsUsbNet^ mcsusb )
11.130.2 Member Function Documentation
11.130.2.1 GetFilterParametersHeadstage() array<int> ^ GetFilterParametersHeadstage (
              unsigned short index )
11.130.2.2 GetHasChecksum() unsigned int GetHasChecksum (
              unsigned int Device )
\textbf{11.130.2.3} \quad \textbf{GetHasRedLedHeadstage()} \quad \texttt{bool GetHasRedLedHeadstage ()} \\
11.130.2.4 GetHeadstageOnOff() USHORT GetHeadstageOnOff ()
\textbf{11.130.2.5} \quad \textbf{GetResetFilter()} \quad \texttt{unsigned int GetResetFilter ()}
              unsigned int Device )
11.130.2.6 GetRFConnectionStatus() unsigned int GetRFConnectionStatus ( )
11.130.2.7 GetRFFrequencyHeadstage() unsigned short GetRFFrequencyHeadstage (
              uint8_t receiver_nb )
```


```
11.130.2.8 GetRFFrequencyReceiver() unsigned short GetRFFrequencyReceiver (
            uint8_t receiver_nb,
            uint8_t configuration )
11.130.2.9 GetRFPower() unsigned short GetRFPower ( )
11.130.2.10 GetScanHeadstagesResult() uint8_t GetScanHeadstagesResult (
            int max_wait_for_ms )
11.130.2.11 GetSelectedHeadstage() uint8_t GetSelectedHeadstage ( )
11.130.2.12 GetSerialNumberHeadstage() unsigned short GetSerialNumberHeadstage ( )
11.130.2.13 GetWPADebugMode() unsigned int GetWPADebugMode (
            unsigned int Device )
11.130.2.14 GetWPAType() unsigned short GetWPAType (
            unsigned int Device )
11.130.2.15 ResetChannelmap() uint32_t ResetChannelmap (
            unsigned int virtualDevice )
11.130.2.16 ScanForHeadstages() void ScanForHeadstages ()
11.130.2.17 SetChannelmap() uint32_t SetChannelmap (
            unsigned char position,
            unsigned char channel,
            unsigned int Device )
```

```
11.130.2.18 SetFilterParametersHeadstage() void SetFilterParametersHeadstage (
              unsigned short index,
              array< int >^{\land} buffer )
11.130.2.19 SetHasChecksum() void SetHasChecksum (
              unsigned int has,
              unsigned int Device )
11.130.2.20 SetHeadstageOnOff() void SetHeadstageOnOff (
              uint16_t onoff )
11.130.2.21 SetHWSelectedChannels() void SetHWSelectedChannels (
              array< bool >^{\wedge} channels,
              unsigned int Device )
11.130.2.22 SetResetFilter() void SetResetFilter (
              unsigned int reset,
              unsigned int Device )
\textbf{11.130.2.23} \quad \textbf{SetRFF} \textbf{requencyHeadstage()} \quad \textbf{void SetRFF} \textbf{requencyHeadstage ()}
              uint8_t receiver_nb,
              unsigned short frequency )
11.130.2.24 SetRFFrequencyReceiver() void SetRFFrequencyReceiver (
              uint8_t receiver_nb,
              uint8_t configuration,
              unsigned short frequency )
11.130.2.25 SetRFFrequencyReceiverEeprom() void SetRFFrequencyReceiverEeprom (
              uint8_t receiver_nb,
              uint8_t configuration,
              unsigned short frequency )
```

```
11.130.2.26 SetRFLostBehaviour() void SetRFLostBehaviour (
            uint8_t stoponfailure,
             unsigned int Device )
11.130.2.27 SetRFPower() void SetRFPower (
             unsigned short power )
11.130.2.28 SetSelectedHeadstage() void SetSelectedHeadstage (
            uint8_t number )
11.130.2.29 SetSerialNumberHeadstage() void SetSerialNumberHeadstage (
            unsigned short number )
11.130.2.30 SetWPADebugMode() void SetWPADebugMode (
             unsigned int mode,
             unsigned int Device )
11.130.2.31 SetWPAType() void SetWPAType (
             unsigned short type,
             unsigned int Device )
```

11.131 CWirelessBaseFunctionNet Class Reference

Inheritance diagram for CWirelessBaseFunctionNet:

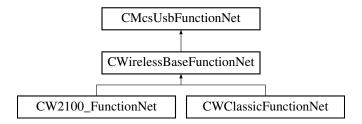

Public Member Functions

• CWirelessBaseFunctionNet (CMcsUsbNet^ mcsusb, CMcsUsbFunctionPointerContainer^ mcsusbfunction)

Static Public Member Functions

• static String ^ CreateWirelessHeadstageSerialNumberString (unsigned short ID)

Additional Inherited Members

11.131.1 Constructor & Destructor Documentation

```
11.131.1.1 CWirelessBaseFunctionNet() CWirelessBaseFunctionNet (
CMcsUsbNet^ mcsusb,
CMcsUsbFunctionPointerContainer^ mcsusbfunction)
```

11.131.2 Member Function Documentation

11.132 DeviceIdNet Struct Reference

Device Id.

Public Member Functions

- DeviceIdNet ()
- DeviceIdNet (VendorIdEnumNet vendor, ProductIdEnumNet product, int bcd, McsBusTypeEnumNet bustype)
- DeviceIdNet (DeviceIdNet% deviceId)
- DeviceIdNet operator= (DeviceIdNet% deviceId)

Public Attributes

- VendorldEnumNet IdVendor
- ProductIdEnumNet IdProduct
- int BcdDevice
- McsBusTypeEnumNet BusType

11.132.1 Detailed Description

Device Id.

11.132.2 Constructor & Destructor Documentation

11.132.3 Member Function Documentation

11.132.2.3 DeviceIdNet() [3/3] DeviceIdNet (

DeviceIdNet% deviceId)

11.132.4 Member Data Documentation

11.132.4.1 BcdDevice int BcdDevice

11.132.4.2 BusType McsBusTypeEnumNet BusType

11.132.4.3 IdProduct ProductIdEnumNet IdProduct

11.132.4.4 IdVendor VendorIdEnumNet IdVendor

11.133 DigitalSource< digitalsourceenum > Class Template Reference

Public Member Functions

- DigitalSource ()
- DigitalSource (digitalsourceenum source)
- int MaxBitNumber ()
- int MaxBitNumber (digitalsourceenum Source)

Static Public Member Functions

- static int MaxBitNumberStatic (digitalsourceenum Source)
- static int size ()

Properties

• digitalsourceenum Source [get, set]

11.133.1 Constructor & Destructor Documentation

```
11.133.1.1 DigitalSource() [1/2] DigitalSource ( )
```

```
11.133.1.2 DigitalSource() [2/2] DigitalSource (
digitalsourceenum source)
```

11.133.2 Member Function Documentation

```
11.133.2.1 MaxBitNumber() [1/2] int MaxBitNumber ( )
```

```
11.133.2.2 MaxBitNumber() [2/2] int MaxBitNumber (
digitalsourceenum Source)
```

```
11.133.2.3 MaxBitNumberStatic() static int MaxBitNumberStatic (
digitalsourceenum Source) [static]
```

```
11.133.2.4 size() static int size () [static]
```

11.133.3 Property Documentation

```
11.133.3.1 Source digitalsourceenum Source [get], [set]
```

11.134 DigitalSourceGeneral Class Reference

Public Member Functions

- DigitalSourceGeneral (Type[^] type)
- DigitalSourceGeneral (Type[^] type, int Source)
- int MaxBitNumber ()
- int MaxBitNumber (int Source)

Static Public Member Functions

- static int MaxBitNumber (Type[^] type, int Source)
- static int size (Type[^] type)

Properties

```
• int Source [get, set]
```

11.134.1 Constructor & Destructor Documentation

```
11.134.1.1 DigitalSourceGeneral() [1/2] DigitalSourceGeneral ( Type^{ \wedge} type )
```

```
11.134.1.2 DigitalSourceGeneral() [2/2] DigitalSourceGeneral ( Type^{\wedge} type, int Source )
```

11.134.2 Member Function Documentation

```
11.134.2.1 MaxBitNumber() [1/3] int MaxBitNumber ( )
11.134.2.2 MaxBitNumber() [2/3] int MaxBitNumber ( int Source )

11.134.2.3 MaxBitNumber() [3/3] static int MaxBitNumber ( Type^ type, int Source ) [static]

11.134.2.4 size() static int size ( Type^ type ) [static]

11.134.3 Property Documentation
```

11.134.3.1 Source int Source [get], [set]

11.135 DriverVersionNet Class Reference

Class gives firmware versions of the device's firmware destinations.

Public Member Functions

• DriverVersionNet ()

Contructor.

∼DriverVersionNet ()

Destructor

unsigned int GetStatus (CFirmwareDestinationNet dest)

Get status of firmware destination.

• unsigned int GetStatus (unsigned int index)

Get status of firmware destination.

unsigned int GetVersionInt (CFirmwareDestinationNet dest)

Get the version number of firmware destination (major in high word, minor in low word)

unsigned int GetVersionInt (unsigned int index)

Get the version number of firmware destination (major in high word, minor in low word)

unsigned int GetMajor (CFirmwareDestinationNet dest)

Get the major version number of firmware destination.

unsigned int GetMajor (unsigned int index)

Get the major version number of firmware destination.

unsigned int GetMinor (CFirmwareDestinationNet dest)

Get the minor version number of firmware destination.

unsigned int GetMinor (unsigned int index)

Get the minor version number of firmware destination.

• unsigned int GetNumEntries ()

Get the number of available firmware destinations.

String \(^\) GetVersionString (CFirmwareDestinationNet dest)

Get the version as a string in the format Major. Minor.

String \(^\) GetVersionString (unsigned int index)

Get the version as a string in the format Major. Minor.

• CFirmwareDestinationNet GetDestinationCode (unsigned int index)

Get CFirmwareDestinationNet.

String \(^\) GetDestinationName (CFirmwareDestinationNet dest)

Get firmware destination name.

String \(^\) GetDestinationName (unsigned int index)

Get firmware destination name.

String \(^\) GetSerialNumber (CFirmwareDestinationNet dest)

Get the serial number of the destination, when no serial number if found, return an empty string.

String \(^\) GetSerialNumber (unsigned int index)

Get the serial number of the destination, when no serial number if found, return an empty string.

Static Public Member Functions

static String ^ DriverVersionNet::FormatVersion (unsigned int v)

11.135.1 Detailed Description

Class gives firmware versions of the device's firmware destinations.

11.135.2 Constructor & Destructor Documentation

```
11.135.2.1 DriverVersionNet() DriverVersionNet ( )
```

Contructor.

11.135.2.2 ~DriverVersionNet() ~DriverVersionNet ()

Destructor.

11.135.3 Member Function Documentation

```
11.135.3.1 DriverVersionNet::FormatVersion() static String ^{\land} DriverVersionNet::FormatVersion ( unsigned int v ) [static]
```

```
11.135.3.2 GetDestinationCode() CFirmwareDestinationNet GetDestinationCode ( unsigned int index )
```

Get CFirmwareDestinationNet.

Parameters

index by index of firmware destination

11.135.3.3 GetDestinationName() [1/2] String $^{\land}$ GetDestinationName (CFirmwareDestinationNet dest)

Get firmware destination name.

Parameters

dest by CFirmwareDestionationNet

11.135.3.4 GetDestinationName() [2/2] String $^{\land}$ GetDestinationName (unsigned int index)

Get firmware destination name.

Parameters

index by index of firmware destination

11.135.3.5 **GetMajor()** [1/2] unsigned int GetMajor (CFirmwareDestinationNet dest)

Get the major version number of firmware destination.

Parameters

dest by CFirmwareDestionationNet

11.135.3.6 **GetMajor()** [2/2] unsigned int GetMajor (unsigned int *index*)

Get the major version number of firmware destination.

Parameters

index	by index of firmware destination

```
11.135.3.7 GetMinor() [1/2] unsigned int GetMinor (
CFirmwareDestinationNet dest)
```

Get the minor version number of firmware destination.

Parameters

```
dest by CFirmwareDestionationNet
```

```
11.135.3.8 GetMinor() [2/2] unsigned int GetMinor ( unsigned int index )
```

Get the minor version number of firmware destination.

Parameters

index by index of firmware destination	ation
--	-------

11.135.3.9 GetNumEntries() unsigned int GetNumEntries ()

Get the number of available firmware destinations.

```
11.135.3.10 GetSerialNumber() [1/2] String ^{\land} GetSerialNumber ( CFirmwareDestinationNet dest )
```

Get the serial number of the destination, when no serial number if found, return an empty string.

Parameters

```
dest by CFirmwareDestionationNet
```

```
11.135.3.11 GetSerialNumber() [2/2] String ^{\land} GetSerialNumber ( unsigned int index )
```

Get the serial number of the destination, when no serial number if found, return an empty string.

Parameters

index by index of firmware destination
--

11.135.3.12 **GetStatus()** [1/2] unsigned int GetStatus (CFirmwareDestinationNet dest)

Get status of firmware destination.

Parameters

11.135.3.13 **GetStatus()** [2/2] unsigned int GetStatus (unsigned int *index*)

Get status of firmware destination.

Parameters

```
11.135.3.14 GetVersionInt() [1/2] unsigned int GetVersionInt ( CFirmwareDestinationNet dest )
```

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

dest	by CFirmwareDestionationNet
------	-----------------------------

```
11.135.3.15 GetVersionInt() [2/2] unsigned int GetVersionInt ( unsigned int index )
```

Get the version number of firmware destination (major in high word, minor in low word)

Parameters

index	by index of firmware destination

```
11.135.3.16 GetVersionString() [1/2] String ^ GetVersionString (
CFirmwareDestinationNet dest)
```

Get the version as a string in the format Major. Minor.

Parameters

```
dest by CFirmwareDestionationNet
```

```
11.135.3.17 GetVersionString() [2/2] String ^{\land} GetVersionString ( unsigned int index )
```

Get the version as a string in the format Major. Minor.

Parameters

index by index of firmware

11.136 FirmwareDestinationNames Class Reference

Static Public Attributes

```
    static String \(^\text{DSP} = \text{gcnew String("DSP")}\)

    static String \(^{\text{USB}} = \text{gcnew String("USB")}\)

    static String ^ MCU1 = gcnew String( "MCU1" )

    static String \(^\text{Bootstrap} = \text{gcnew String("Bootstrap")}\)

    static String \(^{\text{MCSBUS1}} = \text{gcnew String( "McsBus1" )}\)

    static String \(^{\text{MCSBUS2}} = \text{gcnew String( "McsBus2" )}\)

    static String \(^{\text{MCSBUS3}} = \text{gcnew String( "McsBus3" )}\)

    static String \(^\text{MCSBUS4} = \text{gcnew String( "McsBus4" )}\)

    static String \(^{\text{MCSBUS5}} = \text{gcnew String( "McsBus5" )}\)

• static String ^ MCSBUS6 = gcnew String( "McsBus6" )

    static String \(^{\text{MCSBUS7}} = \text{gcnew String( "McsBus7" )}\)

    static String \(^{\text{MCSBUS8}} = \text{gcnew String( "McsBus8" )}\)

    static String \(^{\text{MCSBUS9}} = \text{gcnew String( "McsBus9" )}\)

    static String \(^{\text{MCSBUS10}} = \text{gcnew String( "McsBus10" )}\)

    static String \(^\text{MCSBUS11} = \text{gcnew String("McsBus11")}\)

• static String ^ MCSBUS12 = gcnew String( "McsBus12" )

    static String \(^{\text{MCSBUS13}}\) = gcnew String( "McsBus13" )

    static String \(^\text{BUS1_MCSBUS1} = \text{gcnew String("Bus1McsBus1")}\)

    static String \(^\text{BUS1_MCSBUS2} = \text{gcnew String("Bus1McsBus2")}\)

    static String \(^\text{PIC} = \text{gcnew String("PIC")}\)

    static String \(^\text{PIC2} = \text{gcnew String("PIC2")}\)

    static String \(^\text{PIC3} = \text{gcnew String("PIC3")}\)

• static String ^{\wedge} PIC4 = gcnew String( "PIC4" )
```

```
    static String ^ Altera = gcnew String( "Altera")
    static String ^ FPGA2 = gcnew String( "FPGA2")
    static String ^ FPGA3 = gcnew String( "FPGA3")
    static String ^ FPGA4 = gcnew String( "FPGA4")
    static String ^ FPGA5 = gcnew String( "FPGA5")
    static String ^ FPGA6 = gcnew String( "FPGA6")
```

11.136.1 Member Data Documentation


```
11.136.1.1 Altera String ^ Altera = gcnew String( "Altera" ) [static]
11.136.1.2 Bootstrap String ^ Bootstrap = gcnew String( "Bootstrap") [static]
11.136.1.3 BUS1_MCSBUS1 String ^ BUS1_MCSBUS1 = gcnew String( "BuslMcsBus1" ) [static]
11.136.1.4 BUS1 MCSBUS2 String ^ BUS1_MCSBUS2 = gcnew String( "Bus1McsBus2" ) [static]
11.136.1.5 DSP String ^ DSP = gcnew String( "DSP" ) [static]
11.136.1.6 FPGA2 String ^ FPGA2 = gcnew String( "FPGA2" ) [static]
11.136.1.7 FPGA3 String ^ FPGA3 = gcnew String( "FPGA3" ) [static]
11.136.1.8 FPGA4 String ^ FPGA4 = gcnew String( "FPGA4" ) [static]
11.136.1.9 FPGA5 String ^ FPGA5 = gcnew String( "FPGA5") [static]
```

```
11.136.1.10 FPGA6 String ^ FPGA6 = gcnew String( "FPGA6" ) [static]
11.136.1.11 MCSBUS1 String ^{\wedge} MCSBUS1 = gcnew String( "McsBus1" ) [static]
11.136.1.12 MCSBUS10 String ^ MCSBUS10 = gcnew String( "McsBus10" ) [static]
11.136.1.13 MCSBUS11 String ^ MCSBUS11 = gcnew String( "McsBus11" ) [static]
11.136.1.14 MCSBUS12 String ^ MCSBUS12 = gcnew String( "McsBus12" ) [static]
11.136.1.15 MCSBUS13 String ^ MCSBUS13 = gcnew String( "McsBus13" ) [static]
11.136.1.16 MCSBUS2 String ^ MCSBUS2 = gcnew String( "McsBus2" ) [static]
11.136.1.17 MCSBUS3 String ^{\land} MCSBUS3 = gcnew String( "McsBus3" ) [static]
11.136.1.18 MCSBUS4 String ^ MCSBUS4 = gcnew String( "McsBus4" ) [static]
11.136.1.19 MCSBUS5 String ^ MCSBUS5 = gcnew String( "McsBus5" ) [static]
11.136.1.20 MCSBUS6 String ^ MCSBUS6 = gcnew String( "McsBus6" ) [static]
```

```
11.136.1.21 MCSBUS7 String ^ MCSBUS7 = gcnew String( "McsBus7" ) [static]
11.136.1.22 MCSBUS8 String ^ MCSBUS8 = gcnew String( "McsBus8" ) [static]
11.136.1.23 MCSBUS9 String ^ MCSBUS9 = gcnew String( "McsBus9" ) [static]
11.136.1.24 MCU1 String ^ MCU1 = gcnew String( "MCU1" ) [static]
11.136.1.25 PIC String ^{\wedge} PIC = gcnew String( "PIC" ) [static]
11.136.1.26 PIC2 String ^ PIC2 = gcnew String( "PIC2" ) [static]
11.136.1.27 PIC3 String ^ PIC3 = gcnew String( "PIC3" ) [static]
11.136.1.28 PIC4 String ^ PIC4 = gcnew String( "PIC4" ) [static]
11.136.1.29 USB String ^{\land} USB = gcnew String( "USB" ) [static]
```

11.137 HeadStageIDType Class Reference

Inheritance diagram for HeadStageIDType:

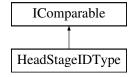

Public Types

enum class HeadstageTypeEnum {
 Unknown ,
 MeasuringOnly ,
 OpticalStimulation ,
 ElectricalStimulation }

Public Member Functions

- HeadStageIDType (unsigned int entry, CW2100 FunctionNet[^] device)
- virtual System::String ^ ToString () override
- virtual bool Equals (Object[^] obj) override
- virtual Int32 CompareTo (Object[^] obj)

Properties

- bool Valid [get]
- unsigned int Entry [get]
- unsigned short ID [get]
- System::String SN [get]
- unsigned int TypeValue [get]
- System::String Type [get]
- HeadstageTypeEnum HeadstageType [get]
- System::String UserDefinedName [get]
- int NumberOfAnalogChannels [get]
- $\bullet \ \, \textbf{int NumberOfStimulationChannels} \quad [\texttt{get}]$
- W2100_StimulusParametersNet^ StimulusParameters [get]
- bool HasIMU [get]
- bool W16lsW14 [get]
- bool HasOptoCurrentMessurement [get]

11.137.1 Member Enumeration Documentation

11.137.1.1 HeadstageTypeEnum enum HeadstageTypeEnum [strong]

Enumerator

Unknown	
MeasuringOnly	
OpticalStimulation	
ElectricalStimulation	

11.137.2 Constructor & Destructor Documentation

```
11.137.2.1 HeadStageIDType() HeadStageIDType (
            unsigned int entry,
            CW2100_FunctionNet^ device )
11.137.3 Member Function Documentation
11.137.3.1 CompareTo() virtual Int32 CompareTo (
            Object^ obj ) [virtual]
11.137.3.2 Equals() virtual bool Equals (
            Object^ obj ) [override], [virtual]
11.137.3.3 ToString() virtual System::String ^ ToString ( ) [override], [virtual]
11.137.4 Property Documentation
11.137.4.1 Entry unsigned int Entry [get]
11.137.4.2 HasIMU bool HasIMU [get]
11.137.4.3 HasOptoCurrentMessurement bool HasOptoCurrentMessurement [get]
11.137.4.4 HeadstageType HeadstageTypeEnum HeadstageType [get]
11.137.4.5 ID unsigned short ID [get]
```

```
11.137.4.6 NumberOfAnalogChannels int NumberOfAnalogChannels [get]
11.137.4.7 NumberOfStimulationChannels int NumberOfStimulationChannels [get]
11.137.4.8 SN System:: String SN [get]
11.137.4.9 StimulusParameters W2100_StimulusParametersNet^ StimulusParameters [get]
11.137.4.10 Type System:: String^{\wedge} Type [get]
11.137.4.11 TypeValue unsigned int TypeValue [get]
11.137.4.12 UserDefinedName System:: String^ UserDefinedName [get]
11.137.4.13 Valid bool Valid [get]
11.137.4.14 W16IsW14 bool W16IsW14 [get]
11.138 HeadstageIDTypeObject Class Reference
```

Public Member Functions

- HeadstageIDTypeObject (HeadStageIDType^ idType)
- virtual String ^ ToString () override
- virtual bool Equals (Object $^{\wedge}$ obj) override
- virtual int GetHashCode () override

Public Attributes

- HeadStageIDType ^ _IdType
- String ^ _AdditionalText

Properties

```
• HeadStageIDType^ IdType [get]
```

```
• String^ AdditionalText [get, set]
```

11.138.1 Constructor & Destructor Documentation

```
11.138.1.1 HeadstageIDTypeObject() HeadstageIDTypeObject ( HeadStageIDType^ idType )
```

11.138.2 Member Function Documentation

```
11.138.2.1 Equals() virtual bool Equals (
          Object^ obj ) [override], [virtual]
```

```
11.138.2.2 GetHashCode() virtual int GetHashCode ( ) [override], [virtual]
```

```
11.138.2.3 ToString() virtual String ^{\wedge} ToString ( ) [override], [virtual]
```

11.138.3 Member Data Documentation

```
\textbf{11.138.3.1} \quad \underline{\textbf{AdditionalText}} \quad \texttt{String} \; \wedge \; \underline{\textbf{AdditionalText}}
```

```
11.138.3.2 _ldType HeadStageIDType ^ _IdType
```

11.138.4 Property Documentation

```
\textbf{11.138.4.1} \quad \textbf{AdditionalText} \quad \texttt{String}^{\wedge} \; \texttt{AdditionalText} \quad \texttt{[get], [set]}
```

```
11.138.4.2 IdType HeadStageIDType^ IdType [get]
```

11.139 HeadStageIDTypeState Class Reference

Properties

- unsigned int State [get]
- HeadStageIDType^ IdType [get]
- bool ControlState [get]
- bool DataState [get]

11.139.1 Property Documentation

```
11.139.1.1 ControlState bool ControlState [get]
```

11.139.1.2 DataState bool DataState [get]

11.139.1.3 IdType HeadStageIDType^ IdType [get]

11.139.1.4 State unsigned int State [get]

11.140 mkfilterNet Class Reference

Static Public Member Functions

- static int mkfilter (String^ filtertype, double value, String^ passtype, int order, double alpha1, double alpha2, [System::Runtime::InteropServices::Out] array< double >^% xcoeffs, [System::Runtime::InteropServices ::Out] array< double >^% ycoeffs)
- static int mkfilter_MCS (int SamplesPerSecond, double R1, double R2, double C, double Amplification, double Correction, [System::Runtime::InteropServices::Out] array< double >^% xcoeffs, [System::Runtime::
 InteropServices::Out] array< double >^% ycoeffs)
- static int mkfilter_MCS (int SamplesPerSecond, double R1, double R2, double C, double Correction, [System::Runtime::InteropServices::Out] array< double >^% xcoeffs, [System::Runtime::InteropServices::Out] array< double >^% ycoeffs)
- static int mkfilter_MCS_k (int SamplesPerSecond, double R1, double R2, double C, double Amplification, double Correction, [System::Runtime::InteropServices::Out] array< double >^% coeffs)
- static int mkfilter_MCS_k (int SamplesPerSecond, double R1, double R2, double C, double Correction, [System::Runtime::InteropServices::Out] array< double >^% coeffs)
- static void mkfilter_coef_in_one_set (int n, [System::Runtime::InteropServices::In] array< double >^ xcoeffs, [System::Runtime::InteropServices::Out] array< double >^% out_coeffs)
- static void mkfilter_scale_coef_in_one_set (int n, double scale, [System::Runtime::InteropServices::In] array< double >^ xcoeffs, [System::Runtime::InteropServices::In] array< double >^ ycoeffs, [System::

 Runtime::InteropServices::Out] array< double >^% out_coeffs)
- static void mkfilter_normalize_coeffs_short (short maxvalue, [System::Runtime::InteropServices::In] array
 double >^ coeffs, [System::Runtime::InteropServices::Out] array< short >^% out_coeffs)
- static void mkfilter_normalize_coeffs_int (int maxvalue, [System::Runtime::InteropServices::In] array< double
 ^ coeffs, [System::Runtime::InteropServices::Out] array< int >^% out_coeffs)
- static void mkfilter_normalize_scale_coeffs_int (int maxvalue, [System::Runtime::InteropServices::In] array
 double >^ coeffs, [System::Runtime::InteropServices::Out] array< int >^% out_coeffs)
- static double mkfilter_highpass_coeff (int SamplesPerSecond, double Frequency)
- static double mkfilter_highpass_k (int SamplesPerSecond, double Frequency)
- static double mkfilter_highpass_frequency_from_coeff (int SamplesPerSecond, double coeff)
- static double mkfilter_highpass_frequency_from_k (int SamplesPerSecond, double k)

11.140.1 Member Function Documentation

```
11.140.1.2 mkfilter_coef_in_one_set() static void mkfilter_coef_in_one_set (
             int n_{i}
             [System::Runtime::InteropServices::In] array< double >^{\wedge} xcoeffs,
             [System::Runtime::InteropServices::In] array< double >^{\wedge} ycoeffs,
             [System::Runtime::InteropServices::Out] array< double >^{\%} out_coeffs ) [static]
11.140.1.3 mkfilter_highpass_coeff() static double mkfilter_highpass_coeff (
             int SamplesPerSecond,
             double Frequency ) [static]
11.140.1.4 mkfilter_highpass_frequency_from_coeff() static double mkfilter_highpass_frequency_←
from coeff (
             int SamplesPerSecond,
             double coeff ) [static]
11.140.1.5 mkfilter highpass frequency from k() static double mkfilter_highpass_frequency_from ↔
_k (
             int SamplesPerSecond,
             double k ) [static]
11.140.1.6 mkfilter_highpass_k() static double mkfilter_highpass_k (
             int SamplesPerSecond,
             double Frequency ) [static]
11.140.1.7 mkfilter_MCS() [1/2] static int mkfilter_MCS (
             int SamplesPerSecond,
             double R1,
             double R2,
             double C,
             double Amplification,
             double Correction,
             [System::Runtime::InteropServices::Out] array< double >^{\%} xcoeffs,
             [System::Runtime::InteropServices::Out] array< double >^% ycoeffs ) [static]
11.140.1.8 mkfilter_MCS() [2/2] static int mkfilter_MCS (
             int SamplesPerSecond,
             double R1,
             double R2,
             double C_{\prime}
             double Correction,
             [System::Runtime::InteropServices::Out] array< double >^{\%} xcoeffs,
             [System::Runtime::InteropServices::Out] array< double >^{\%} ycoeffs ) [static]
```

```
11.140.1.9 mkfilter_MCS_k() [1/2] static int mkfilter_MCS_k (
              int SamplesPerSecond,
              double R1,
              double R2,
              double C,
              double Amplification,
              double Correction,
              [System::Runtime::InteropServices::Out] array< double >^{\%} coeffs ) [static]
11.140.1.10 mkfilter_MCS_k() [2/2] static int mkfilter_MCS_k (
              int SamplesPerSecond,
              double R1.
              double R2,
              double C,
              double Correction,
              [System::Runtime::InteropServices::Out] array< double >^{\%} coeffs ) [static]
\textbf{11.140.1.11} \quad \textbf{mkfilter\_normalize\_coeffs\_int()} \quad \texttt{static void mkfilter\_normalize\_coeffs\_int} \ (
              int maxvalue,
              [System::Runtime::InteropServices::In] array< double >^{\wedge} coeffs,
              [System::Runtime::InteropServices::Out] array< int >^{\%} out_coeffs ) [static]
11.140.1.12 mkfilter_normalize_coeffs_short() static void mkfilter_normalize_coeffs_short (
              short maxvalue,
              [System::Runtime::InteropServices::In] array< double >^{\wedge} coeffs,
              [System::Runtime::InteropServices::Out] array< short >^{\%} out_coeffs ) [static]
11.140.1.13 mkfilter_normalize_scale_coeffs_int() static void mkfilter_normalize_scale_coeffs_int
              int maxvalue,
              [System::Runtime::InteropServices::In] array< double >^{\land} coeffs,
               [System::Runtime::InteropServices::Out] \ array< int > ^\$ \ out\_coeffs \ ) \ [static] 
11.140.1.14 mkfilter_scale_coef_in_one_set() static void mkfilter_scale_coef_in_one_set (
              int n,
              double scale,
              [System::Runtime::InteropServices::In] array< double >^{\wedge} xcoeffs,
              [System::Runtime::InteropServices::In] array< double >^{\wedge} ycoeffs,
              [System::Runtime::InteropServices::Out] \ array< \ double > ^ {\$} \ out\_coeffs \ ) \ [static]
```

11.141 CRoboDeviceNet::RoboMainLowLevelCommands Class Reference

Public Member Functions

- void SetParameter (unsigned short command, unsigned short index, unsigned int value)
- void SetParameter (unsigned short command, unsigned short index, unsigned int value1, unsigned int value2)
- void SetUserParameter (unsigned short index, unsigned int value)

Stores persistently 32 bit integer values on RoboMain

void SetUserParameter (unsigned short index, int value)

Stores persistently 32 bit integer values on RoboMain

- void GetParameter (unsigned short command, unsigned short index, [System::Runtime::InteropServices::
 — Out]unsigned int% value)
- void GetParameter (unsigned short command, unsigned short index, [System::Runtime::InteropServices::
 Out]unsigned int% value1, [System::Runtime::InteropServices::Out]unsigned int% value2)
- void GetUserParameter (unsigned short index, [System::Runtime::InteropServices::Out]unsigned int% value)

Reads 32 bit integer values stored persistently on RoboMain

void GetUserParameter (unsigned short index, [System::Runtime::InteropServices::Out]int% value)

Reads 32 bit integer values stored persistently on RoboMain

- void FindReferencePhase0 (unsigned char busaddress, char axes)
- · void FindReferencePhase0 (unsigned char busaddress, char axes, int timeout)
- unsigned char HasRef (unsigned char busaddress, char axes)
- void SetHWRevision (unsigned int revision)
- unsigned int GetHWRevision ()
- · void SetHWConfig (unsigned int config)
- unsigned int GetHWConfig ()
- void SetMinPressureWaitTime (unsigned int t)
- unsigned int GetMinPressureWaitTime ()
- void SetMinPressure (unsigned int pressure)
- unsigned int GetMinPressure ()
- void SetMaxPressureWaitTime (unsigned int t)
- unsigned int GetMaxPressureWaitTime ()
- void SetMinNoPressureWaitTime (unsigned int t)
- unsigned int GetMinNoPressureWaitTime ()
- void SetMaxNoPressure (unsigned int pressure)
- unsigned int GetMaxNoPressure ()
- void SetMaxNoPressureWaitTime (unsigned int t)
- unsigned int GetMaxNoPressureWaitTime ()
- · void SetSearchReferenceMethod (unsigned char busaddress, char axes, unsigned int method)
- unsigned int GetSearchReferenceMethod (unsigned char busaddress, char axes)
- void SetSearchReferenceOffsetPos (unsigned char busaddress, char axes, int offsetpos)
- int GetSearchReferenceOffsetPos (unsigned char busaddress, char axes)
- void SetSearchReferenceFastSpeed (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short GetSearchReferenceFastSpeed (unsigned char busaddress, char axes)
- · void SetSearchReferenceFastAccel (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short GetSearchReferenceFastAccel (unsigned char busaddress, char axes)
- void SetSearchReferenceFineSpeed (unsigned char busaddress, char axes, unsigned short speed)
- unsigned short GetSearchReferenceFineSpeed (unsigned char busaddress, char axes)
- · void SetSearchReferenceFineAccel (unsigned char busaddress, char axes, unsigned short accel)
- unsigned short GetSearchReferenceFineAccel (unsigned char busaddress, char axes)
- void SetSearchReferenceMoveOut (unsigned char busaddress, char axes, int move)
- int GetSearchReferenceMoveOut (unsigned char busaddress, char axes)
- void SetAxisConfig (unsigned char busaddress, char axes, unsigned int config)
- unsigned int GetAxisConfig (unsigned char busaddress, char axes)
- void GetPhases (unsigned char busaddress, char axes, [System::Runtime::InteropServices::Out] unsigned short% phase0, [System::Runtime::InteropServices::Out] unsigned short% lastphase)

11.141.1 Member Function Documentation

```
11.141.1.1 FindReferencePhaseO() [1/2] void FindReferencePhaseO (
             unsigned char busaddress,
             char axes )
\textbf{11.141.1.2} \quad \textbf{FindReferencePhase0() [2/2]} \quad \texttt{void FindReferencePhase0} \quad \textbf{(}
             unsigned char busaddress,
             char axes,
             int timeout )
11.141.1.3 GetAxisConfig() unsigned int GetAxisConfig (
             unsigned char busaddress,
             char axes )
11.141.1.4 GetHWConfig() unsigned int GetHWConfig ( )
11.141.1.5 GetHWRevision() unsigned int GetHWRevision ()
11.141.1.6 GetMaxNoPressure() unsigned int GetMaxNoPressure ( )
11.141.1.7 GetMaxNoPressureWaitTime() unsigned int GetMaxNoPressureWaitTime ( )
11.141.1.8 GetMaxPressureWaitTime() unsigned int GetMaxPressureWaitTime ( )
11.141.1.9 GetMinNoPressureWaitTime() unsigned int GetMinNoPressureWaitTime ( )
```

```
11.141.1.10 GetMinPressure() unsigned int GetMinPressure ( )
11.141.1.11 GetMinPressureWaitTime() unsigned int GetMinPressureWaitTime ( )
11.141.1.12 GetParameter() [1/2] void GetParameter (
             unsigned short command,
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned int% value )
11.141.1.13 GetParameter() [2/2] void GetParameter (
             unsigned short command,
             unsigned short index,
             [System::Runtime::InteropServices::Out] unsigned int% value1,
             [System::Runtime::InteropServices::Out] unsigned int% value2 )
11.141.1.14 GetPhases() void GetPhases (
             unsigned char busaddress,
             char axes,
             [System::Runtime::InteropServices::Out] unsigned short% phase0,
             [System::Runtime::InteropServices::Out] unsigned short% lastphase )
11.141.1.15 GetSearchReferenceFastAccel() unsigned short GetSearchReferenceFastAccel (
             unsigned char busaddress,
             char axes )
11.141.1.16 GetSearchReferenceFastSpeed() unsigned short GetSearchReferenceFastSpeed (
             unsigned char busaddress,
             char axes )
\textbf{11.141.1.17} \quad \textbf{GetSearchReferenceFineAccel()} \quad \textbf{unsigned short GetSearchReferenceFineAccel ()}
             unsigned char busaddress,
             char axes )
```

```
11.141.1.18 GetSearchReferenceFineSpeed() unsigned short GetSearchReferenceFineSpeed (unsigned char busaddress, char axes)

11.141.1.19 GetSearchReferenceMethod() unsigned int GetSearchReferenceMethod (unsigned char busaddress, char axes)

11.141.1.20 GetSearchReferenceMoveOut() int GetSearchReferenceMoveOut (unsigned char busaddress, char axes)

11.141.1.21 GetSearchReferenceOffsetPos() unsigned char busaddress, char axes)

11.141.1.22 GetUserParameter() [1/2] void GetUserParameter (unsigned short index,
```

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

index	address offset of parameter; range: 015
value	data buffer

```
11.141.1.23 GetUserParameter() [2/2] void GetUserParameter (
unsigned short index,

[System::Runtime::InteropServices::Out] unsigned int% value )
```

[System::Runtime::InteropServices::Out] int% value)

Reads 32 bit integer values stored persistently on RoboMain

intention: provide free persistent user memory space on motor controller

Parameters

index	address offset of parameter; range: 015
value	data buffer

```
11.141.1.24 HasRef() unsigned char HasRef (
            unsigned char busaddress,
             char axes )
11.141.1.25 SetAxisConfig() void SetAxisConfig (
            unsigned char busaddress,
             char axes,
             unsigned int config )
11.141.1.26 SetHWConfig() void SetHWConfig (
            unsigned int config )
11.141.1.27 SetHWRevision() void SetHWRevision (
             unsigned int revision )
11.141.1.28 SetMaxNoPressure() void SetMaxNoPressure (
             unsigned int pressure )
11.141.1.29 SetMaxNoPressureWaitTime() void SetMaxNoPressureWaitTime (
             unsigned int t )
11.141.1.30 SetMaxPressureWaitTime() void SetMaxPressureWaitTime (
            unsigned int t )
11.141.1.31 SetMinNoPressureWaitTime() void SetMinNoPressureWaitTime (
            unsigned int t )
11.141.1.32 SetMinPressure() void SetMinPressure (
             unsigned int pressure )
```

```
11.141.1.33 SetMinPressureWaitTime() void SetMinPressureWaitTime (
             unsigned int t )
11.141.1.34 SetParameter() [1/2] void SetParameter (
             unsigned short command,
             unsigned short index,
             unsigned int value )
11.141.1.35 SetParameter() [2/2] void SetParameter (
             unsigned short command,
             unsigned short index,
             unsigned int value1,
             unsigned int value2 )
11.141.1.36 SetSearchReferenceFastAccel() void SetSearchReferenceFastAccel (
             unsigned char busaddress,
             char axes,
             unsigned short accel )
11.141.1.37 SetSearchReferenceFastSpeed() void SetSearchReferenceFastSpeed (
             unsigned char busaddress,
             char axes,
             unsigned short speed )
11.141.1.38 SetSearchReferenceFineAccel() void SetSearchReferenceFineAccel (
             unsigned char busaddress,
             char axes,
             unsigned short accel )
11.141.1.39 SetSearchReferenceFineSpeed() void SetSearchReferenceFineSpeed (
             unsigned char busaddress,
             char axes,
             unsigned short speed )
```

```
11.141.1.40 SetSearchReferenceMethod() void SetSearchReferenceMethod (
```

```
unsigned char busaddress,
char axes,
unsigned int method )
```

11.141.1.41 SetSearchReferenceMoveOut() void SetSearchReferenceMoveOut (

```
unsigned char busaddress,
char axes,
int move )
```

11.141.1.42 SetSearchReferenceOffsetPos() void SetSearchReferenceOffsetPos (

```
unsigned char busaddress,
char axes,
int offsetpos )
```

11.141.1.43 SetUserParameter() [1/2] void SetUserParameter (

```
unsigned short index,
int value )
```

Stores persistently 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

index	address offset of parameter; range: 015
value	data to be stored

11.141.1.44 SetUserParameter() [2/2] void SetUserParameter (unsigned short index,

```
unsigned short index, unsigned int value )
```

Stores persistently 32 bit integer values on RoboMain

intention: provide free persistent user memory space on RoboMain

Parameters

index	address offset of parameter; range: 015
value	data to be stored

11.142 CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands Class Reference

Public Member Functions

- void FindReferencePhase0XY ()
- void FindReferencePhase0XY (int timeout)

11.142.1 Member Function Documentation

```
11.142.1.1 FindReferencePhase0XY() [1/2] void FindReferencePhase0XY ( )
```

```
11.142.1.2 FindReferencePhaseOXY() [2/2] void FindReferencePhaseOXY ( int timeout )
```

11.143 CFilterCoefficientsNet::s_FilterAttributesNet Struct Reference

Public Member Functions

- s_FilterAttributesNet (s_FilterAttributes attrib)
- s_FilterAttributes ToCpp ()

Public Attributes

- uint32_t PreCommaB
- uint32_t PostCommaB
- uint32_t CommaPositionB
- uint32_t PreCommaA
- uint32_t PostCommaA
- uint32_t CommaPositionA

11.143.1 Constructor & Destructor Documentation

```
11.143.1.1 s_FilterAttributesNet() s_FilterAttributesNet( s_FilterAttributes attrib)
```

11.143.2 Member Function Documentation

11.143.2.1 ToCpp() s_FilterAttributes ToCpp ()
11.143.3 Member Data Documentation
11.143.3.1 CommaPositionA uint32_t CommaPositionA
11.143.3.2 CommaPositionB uint32_t CommaPositionB
11.143.3.3 PostCommaA uint32_t PostCommaA
11.143.3.4 PostCommaB uint32_t PostCommaB
11.143.3.5 PreCommaA uint32_t PreCommaA
11.143.3.6 PreCommaB uint32_t PreCommaB
11.144 CMeaAudioFunctionNet::s_setaudionet Struct Reference
Public Attributes
int channelint amplification
11.144.1 Member Data Documentation
11.144.1.1 amplification int amplification

```
11.144.1.2 channel int channel
```

11.145 CStimulusFunctionNet::SidebandData Class Reference

Public Member Functions

- SidebandData ()
- \sim SidebandData ()

Destructor: called by Dispose()

· !SidebandData ()

Finalizer: called by GC before collecting

Properties

```
• array< int32_t >^ Sideband [get]
```

• array< uint64_t >^ Duration [get]

11.145.1 Constructor & Destructor Documentation

```
11.145.1.1 SidebandData() SidebandData ()
```

```
11.145.1.2 ~SidebandData() ~SidebandData ()
```

Destructor: called by Dispose()

```
11.145.1.3 "!SidebandData() !SidebandData ()
```

Finalizer: called by GC before collecting

11.145.2 Property Documentation

```
11.145.2.1 Duration array< uint64_t>^{\wedge} Duration [get]
```

11.145.2.2 Sideband array< int32_t> $^{\land}$ Sideband [get]

11.146 StgStatusNet Class Reference

Static Public Member Functions

- static StgStatusNet ^ FromIntPtr (IntPtr stgstatus)
- static StgStatusNet ^ FromPtr (stgstatus_t *stgstatus)

Public Attributes

- array< Stg200xTriggerStatusEnumNet > ^ TiggerStatus
- array< uint32_t > ^ ListOfChangedTriggers

11.146.1 Member Function Documentation

```
11.146.1.2 FromPtr() static StgStatusNet ^ FromPtr (
    stgstatus_t * stgstatus ) [static]
```

11.146.2 Member Data Documentation

```
\textbf{11.146.2.1} \quad \textbf{ListOfChangedTriggers} \quad \texttt{array} < \texttt{uint32\_t} > \\ \land \quad \texttt{ListOfChangedTriggers}
```

```
11.146.2.2 TiggerStatus array<Stg200xTriggerStatusEnumNet> ^ TiggerStatus
```

11.147 CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData Class Reference

Public Member Functions

- StimulusDeviceDataAndUnrolledData ()
- ~StimulusDeviceDataAndUnrolledData ()

Destructor: called by Dispose()

• !StimulusDeviceDataAndUnrolledData ()

Finalizer: called by GC before collecting

Properties

• array< uint8_t >^ DeviceData [get]

```
• int DeviceDataLength [get]
   • array< int32_t >^ UnrolledAmplitude [get]
   • array< uint32_t >^ UnrolledSync [get]
   • array< uint64_t >^ UnrolledDuration [get]
11.147.1 Constructor & Destructor Documentation
11.147.1.1 StimulusDeviceDataAndUnrolledData() StimulusDeviceDataAndUnrolledData ()
11.147.1.2 ~StimulusDeviceDataAndUnrolledData() ~StimulusDeviceDataAndUnrolledData ()
Destructor: called by Dispose()
11.147.1.3 "!StimulusDeviceDataAndUnrolledData() !StimulusDeviceDataAndUnrolledData ()
Finalizer: called by GC before collecting
11.147.2 Property Documentation
11.147.2.1 DeviceData array< uint8_t>^ DeviceData [get]
11.147.2.2 DeviceDataLength int DeviceDataLength [get]
11.147.2.3 UnrolledAmplitude array< int32_t>^ UnrolledAmplitude [get]
11.147.2.4 UnrolledDuration array< uint64_t>^ UnrolledDuration [get]
```

11.147.2.5 UnrolledSync array< uint32_t>^ UnrolledSync [get]

11.148 CMcsUsbFactoryNet::TargetSelectionInfoNet Class Reference

Public Attributes

- uint32 t Address
- uint32_t Value

11.148.1 Member Data Documentation

```
11.148.1.1 Address uint32_t Address
```

```
11.148.1.2 Value uint32_t Value
```

11.149 usbSetupPacket_t Class Reference

Public Attributes

- uint8_t bmRequestType
- uint8_t bRequest
- uint16_t wValue
- uint16_t wIndex
- uint16_t wLength

11.149.1 Member Data Documentation

```
11.149.1.1 bmRequestType uint8_t bmRequestType
```

11.149.1.2 bRequest uint8_t bRequest

11.149.1.3 wlndex uint16_t wIndex

- 11.149.1.4 wLength uint16_t wLength
- 11.149.1.5 wValue uint16_t wValue

11.150 W2100_StimulusParametersNet Struct Reference

Public Attributes

- int DACResolution
- int TimeResolutionInNanoSeconds
- int VoltageRangeInMicroVolt
- int VoltageResolutionInMicroVolt
- int CurrentRangeInNanoAmp
- int CurrentResolutionInNanoAmp

11.150.1 Member Data Documentation

- 11.150.1.1 CurrentRangeInNanoAmp int CurrentRangeInNanoAmp
- 11.150.1.2 CurrentResolutionInNanoAmp int CurrentResolutionInNanoAmp
- 11.150.1.3 DACResolution int DACResolution
- $\textbf{11.150.1.4} \quad \textbf{TimeResolutionInNanoSeconds} \quad \texttt{int TimeResolutionInNanoSeconds}$
- $\textbf{11.150.1.5} \quad \textbf{VoltageRangeInMicroVolt} \quad \texttt{int VoltageRangeInMicroVolt}$
- 11.150.1.6 VoltageResolutionInMicroVolt int VoltageResolutionInMicroVolt

Index

!CDacCalibrationFunctionNet	CUsbDeviceConfigurationFunctionNet, 612
CDacCalibrationFunctionNet, 119	!CWarnerUssingDeviceNet
!CDigOutStimulatorFunctionNet	CWarnerUssingDeviceNet, 630
CDigOutStimulatorFunctionNet, 129	!CWarnerUssingFunctionNet
!CExternDTesterDeviceNet	CWarnerUssingFunctionNet, 632
CExternDTesterDeviceNet, 134	!CWarnerValveControllerDeviceNet
!CGrapheneFunctionNet	CWarnerValveControllerDeviceNet, 651
CGrapheneFunctionNet, 173	!CWarnerValveControllerDeviceTesterFunctionNet
!CInterfaceboard2FunctionNet	CWarnerValveControllerDeviceTesterFunctionNet,
CInterfaceboard2FunctionNet, 192	673
!CInterfaceboardFunctionNet	!SidebandData
CInterfaceboardFunctionNet, 194	CStimulusFunctionNet::SidebandData, 711
!CLIH3DeviceNet	!StimulusDeviceDataAndUnrolledData
CLIH3DeviceNet, 197	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
!CMEA2100x256FunctionNet	713
CMEA2100x256FunctionNet, 340	AdditionalText
!CMcsUsbFunctionNet	HeadstageIDTypeObject, 697
CMcsUsbFunctionNet, 309	_ldType
!CMcsUsbListNet	HeadstageIDTypeObject, 697
CMcsUsbListNet, 316	~CCMOSMeaDeviceNet
!CMcsUsbNet	CCMOSMeaDeviceNet, 114
CMcsUsbNet, 322	\sim CChannelTestDeviceNet
!CMeFunctionNet	CChannelTestDeviceNet, 101
CMeFunctionNet, 374	~CCreateFilterNet
!CMeaCleanDeviceNet	CCreateFilterNet, 117
CMeaCleanDeviceNet, 345	~CDacCalibrationFunctionNet
!CMeaCoatDeviceNet	CDacCalibrationFunctionNet, 119
CMeaCoatDeviceNet, 349	~CDigOutStimulatorFunctionNet
!CMultiBatteryChargerDeviceNet	CDigOutStimulatorFunctionNet, 129
CMultiBatteryChargerDeviceNet, 376	~CExternDTesterDeviceNet
!CMultiwellCallbackFunctionNet	CExternDTesterDeviceNet, 134
CMultiwellCallbackFunctionNet, 383	~CFilterCoefficientsNet
!CMultiwellDeviceNet	
CMultiwellDeviceNet, 385	CFilterPrepartyNet
	~CFilterPropertyNet
!CMultiwellOptoStimFunctionNet	CFilterPropertyNet, 141
CMultiwellOptoStimFunctionNet, 392 !CPPCFunctionNet	~CFluidControlDeviceNet
	CFluidControlDeviceNet, 144
CPPCFunctionNet, 426	~CGenericDevelopDeviceNet
!CPedoterDeviceNet	CGenericDevelopDeviceNet, 157
CPedoterDeviceNet, 407	~CGilsonDeviceNet
!CPositionIIDeviceNet	CGilsonDeviceNet, 169
CPositionIIDeviceNet, 413	~CGrapheneASICDeviceNet
!CPositionImpDeviceNet	CGrapheneASICDeviceNet, 170
CPositionImpDeviceNet, 421	~CGrapheneFunctionNet
!CProgramPressureCurveNet	CGrapheneFunctionNet, 173
CProgramPressureCurveNet, 440	~CInterfaceboard2FunctionNet
!CPulseGeneratorFunctionNet	CInterfaceboard2FunctionNet, 192
CPulseGeneratorFunctionNet, 441	~CInterfaceboardFunctionNet
!CRFFunctionNet	CInterfaceboardFunctionNet, 194
CRFFunctionNet, 449	~CLIH3DeviceNet
!CSCUFunctionNet	CLIH3DeviceNet, 197
CSCUFunctionNet, 505	~CMEA2100x256FunctionNet
!CTEERFunctionNet	CMEA2100x256FunctionNet, 340
CTEERFunctionNet, 602	~CMcsBusNet
!CUsbDeviceConfigurationFunctionNet	CMcsBusNet, 243

~CMcsBus AxisParametersNet	~CPeristalticPumpDeviceNet
CMcsBus_AxisParametersNet, 205	CPeristalticPumpDeviceNet, 408
~CMcsBus ExtensionNet	~CPgaDeviceNet
CMcsBus_ExtensionNet, 206	CPgaDeviceNet, 409
~CMcsBus_FYIExtensionNet	\sim CPositionIIDeviceNet
CMcsBus_FYIExtensionNet, 207	CPositionIIDeviceNet, 413
~CMcsBus_MotorControlNet	~CPositionImpDeviceNet
CMcsBus_MotorControlNet, 212	CPositionImpDeviceNet, 421
~CMcsBus SensorNet	~CProgramPressureCurveNet
CMcsBus SensorNet, 228	CProgramPressureCurveNet, 439
~CMcsBus_TempSensorNet	~CPulseGeneratorFunctionNet
CMcsBus TempSensorNet, 237	CPulseGeneratorFunctionNet, 441
~CMcsBus_VoltageModeNet	\sim CRFFunctionNet
CMcsBus_VoltageModeNet, 239	CRFFunctionNet, 449
~CMcsUsbDacqNet	\sim CRetinaLedDeviceNet
CMcsUsbDacqNet, 253	CRetinaLedDeviceNet, 447
~CMcsUsbFactoryNet	~CRoboDeviceNet
CMcsUsbFactoryNet, 301	CRoboDeviceNet, 474
~CMcsUsbFunctionNet	\sim CRoboFluidDeviceNet
CMcsUsbFunctionNet, 309	CRoboFluidDeviceNet, 488
~CMcsUsbListEntryNet	\sim CSCUFunctionNet
CMcsUsbListEntryNet, 310	CSCUFunctionNet, 505
~CMcsUsbListNet	\sim CSafeISDeviceNet
CMcsUsbListNet, 316	CSafeISDeviceNet, 499
~CMcsUsbNet	~CStg200xBasicNet
CMcsUsbNet, 322	CStg200xBasicNet, 525
~CMeFunctionNet	~CStg200xDownloadNet
CMeFunctionNet, 374	CStg200xDownloadNet, 568
~CMeaCleanDeviceNet	~CSw2to64DeviceNet
CMeaCleanDeviceNet, 345	CSw2to64DeviceNet, 585
~CMeaCoatDeviceNet	\sim CTEERFunctionNet
CMeaCoatDeviceNet, 349	CTEERFunctionNet, 602
~CMeaDeviceNet	\sim CTEERMachineDeviceNet
CMeaDeviceNet, 355	CTEERMachineDeviceNet, 610
~CMeaImpedanceDeviceNet	\sim CTcxDeviceNet
CMealmpedanceDeviceNet, 368	CTcxDeviceNet, 589
~CMeaSwitchDeviceNet	\sim CUsbDeviceConfigurationFunctionNet
CMeaSwitchDeviceNet, 371	CUsbDeviceConfigurationFunctionNet, 612
~CMeaUSBDeviceNet	~CWarnerUssingDeviceNet
CMeaUSBDeviceNet, 373	CWarnerUssingDeviceNet, 629
~CMultiBatteryChargerDeviceNet	~CWarnerUssingFunctionNet
CMultiBatteryChargerDeviceNet, 376	CWarnerUssingFunctionNet, 632
~CMultiwellCallbackFunctionNet	\sim CWarnerValveControllerDeviceNet
CMultiwellCallbackFunctionNet, 383	CWarnerValveControllerDeviceNet, 651
~CMultiwellDeviceNet	\sim CWarnerValveControllerDeviceTesterFunctionNet
CMultiwellDeviceNet, 385	CWarnerValveControllerDeviceTesterFunctionNet,
~CMultiwellOptoStimFunctionNet	672
CMultiwellOptoStimFunctionNet, 392	\sim DriverVersionNet
~CNF_GenDeviceNet	DriverVersionNet, 685
CNF_GenDeviceNet, 396	\sim SidebandData
~COkuvisionStimulatorDeviceNet	CStimulusFunctionNet::SidebandData, 711
COkuvisionStimulatorDeviceNet, 402	\sim StimulusDeviceDataAndUnrolledData
\sim CPPCFunctionNet	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
CPPCFunctionNet, 426	713
\sim CPathIdentDeviceNet	
CPathIdentDeviceNet, 406	Α
\sim CPedoterDeviceNet	CFilterCoefficientsNet, 137
CPedoterDeviceNet, 407	Mcs::Usb, 69
	AccelOnly

Mcs::Usb, 91	AnalogSource_IF
AdapterTypeEnumNet	Mcs::Usb, 53
Mcs::Usb, 52	AnalogSourceEnumNet
AdditionalText	Mcs::Usb, 53
HeadstageIDTypeObject, 698	AnalogUnitEnumNet
AddLoopEntry	Mcs::Usb, 53
CRetinaLedDeviceNet, 447	Any
Address	Mcs::Usb, 77, 91
CMcsUsbFactoryNet::TargetSelectionInfoNet, 714	ApplyGains
AddSelectedChannelsQueue	CPgaDeviceNet, 410
CMcsUsbDacqNet, 253-256	AreTransistorVoltagesSet
AddSoftwareKey	CCMOSMea FunctionNet, 104
CMcsUsbNet, 322	Armed
AddTableEntry	Mcs::Usb, 85
CRetinaLedDeviceNet, 447	ASMedia
ALA VC3	Mcs::Usb, 90
Mcs::Usb, 77, 91	AssociateToThis
ALA_VC3_DEVICE	CMcsUsbNet, 322
Mcs::Usb, 63	AudioTestChannelGroup
ALTERA	Mcs::Usb, 61, 71, 82, 93
Mcs::Usb, 55	AutomaticAnalogOut
Altera	CSCUFunctionNet, 505
FirmwareDestinationNames, 691	Aux
ALTERA BASE	Mcs::Usb, 63
_	AuxIn
Mcs::Usb, 56	_
ALTERA_BOOTSTRAP	Mcs::Usb, 57, 71, 73, 74, 83, 86, 93
Mcs::Usb, 56	AuxPort
ALTERA_GOLD	Mcs::Usb, 59
Mcs::Usb, 55	Axes_I
ALTERA_TARGET1	CRoboDeviceNet, 482
Mcs::Usb, 56	Axes_X
ALTERA_TARGET2	CRoboDeviceNet, 482
Mcs::Usb, 56	Axes_Y
ALTERA_TARGET3	CRoboDeviceNet, 482
Mcs::Usb, 56	Axes_Z
AlwaysOn	CRoboDeviceNet, 483
Mcs::Usb, 79	Axis_I
Ampere	CRoboDeviceNet, 483
Mcs::Usb, 53, 91	Axis_X
amplification	CRoboDeviceNet, 483
CMeaAudioFunctionNet::s_setaudionet, 710	Axis_Y
CW2100_FunctionNet::AudioChannelsNet, 96	CRoboDeviceNet, 483
AmplifierSettle	Axis_Z
CIntanMea_FunctionNet, 190	CRoboDeviceNet, 483
AMS_Dongle	D
Mcs::Usb, 77	В
Analog	CFilterCoefficientsNet, 137
Mcs::Usb, 76, 94	Mcs::Usb, 69
AnalogGain	BatteryState, 96
CMeaDeviceNet, 360	Charge, 96
AnalogGroup	ChargeRegionString, 96
Mcs::Usb, 61	ChargeString, 96
AnalogOut_DAC_Range_EnumNet	Voltage, 97
Mcs::Usb, 52	VoltageString, 97
AnalogSource_HS1	BcdDevice
Mcs::Usb, 53	DeviceIdNet, 681
AnalogSource_HS2	BeginImpedanceCheck
Mcs::Usb, 53	CIntanMea_FunctionNet, 190
	Bessel

M	NA 111 54
Mcs::Usb, 66	Mcs::Usb, 54
BesselFilterBandPassNet, 97	BUSOMCSBUS9
BesselFilterBandPassNet, 97	Mcs::Usb, 54
BesselFilterHighPassNet, 97	BUS1_MCSBUS1
BesselFilterHighPassNet, 98	FirmwareDestinationNames, 691
BesselFilterLowPassNet, 98	BUS1_MCSBUS2
BesselFilterLowPassNet, 98 BMI	FirmwareDestinationNames, 691
	BUS1MCSBUS0
Mcs::Usb, 87	Mcs::Usb, 54 BUS1MCSBUS1
bmRequestType	
usbSetupPacket_t, 714 BOOST_BIT	Mcs::Usb, 54 BUS1MCSBUS10
CW2100_StimulatorFunctionNet, 628	Mcs::Usb, 54
	BUS1MCSBUS11
Bootstrap FirmwareDestinationNames, 691	Mcs::Usb, 54
Mcs::Usb, 55	BUS1MCSBUS12
BootstrapOtherCypress	Mcs::Usb, 54
Mcs::Usb, 55	BUS1MCSBUS13
Both	Mcs::Usb, 54
Mcs::Usb, 91	BUS1MCSBUS14
Break	Mcs::Usb, 54
Mcs::Usb, 81	BUS1MCSBUS15
bRequest	Mcs::Usb, 54
usbSetupPacket_t, 714	BUS1MCSBUS2
BurnAdcOffset	Mcs::Usb, 54
COctoPotDeviceNet, 398	BUS1MCSBUS3
BurnDacOffset	Mcs::Usb, 54
CDacCalibrationFunctionNet, 119	BUS1MCSBUS4
COctoPotDeviceNet, 398	Mcs::Usb, 54
BUSOMCSBUSO	BUS1MCSBUS5
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS1	BUS1MCSBUS6
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS10	BUS1MCSBUS7
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS11	BUS1MCSBUS8
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS12	BUS1MCSBUS9
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS13	BUS2MCSBUS0
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS14	BUS2MCSBUS1
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS15	BUS2MCSBUS10
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS2	BUS2MCSBUS11
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS3	BUS2MCSBUS12
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS4	BUS2MCSBUS13
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS5	BUS2MCSBUS14
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS6	BUS2MCSBUS15
Mcs::Usb, 54	Mcs::Usb, 55
BUS0MCSBUS7	BUS2MCSBUS2
Mcs::Usb, 54	Mcs::Usb, 54
BUS0MCSBUS8	BUS2MCSBUS3

Mcs::Usb, 54	CatchAmpGetDacEnable
BUS2MCSBUS4	CMcsBus_SensorNet, 229
Mcs::Usb, 54	CatchAmpGetDacOffset
BUS2MCSBUS5	CMcsBus_SensorNet, 229
Mcs::Usb, 55	CatchAmpGetPwmEnable
BUS2MCSBUS6	CMcsBus_SensorNet, 229
Mcs::Usb, 55	CatchAmpSetDacAmplitude
BUS2MCSBUS7	CMcsBus_SensorNet, 229
Mcs::Usb, 55	CatchAmpSetDacEnable
BUS2MCSBUS8	CMcsBus_SensorNet, 229
Mcs::Usb, 55	CatchAmpSetDacOffset
BUS2MCSBUS9	CMcsBus_SensorNet, 229
Mcs::Usb, 55	CatchAmpSetPwmEnable
BUSNUMBER0	CMcsBus_SensorNet, 229
Mcs::Usb, 54	CChannelTestDeviceNet, 100
BUSNUMBER1	~CChannelTestDeviceNet, 101
Mcs::Usb, 54	CChannelTestDeviceNet, 101
BUSNUMBER2	SetAmplitude, 101
Mcs::Usb, 54	SetAttenuation, 101
BusType	SetFrequency, 101
DeviceIdNet, 681	SetWaveform, 101
Butterworth	CCMOSMea_FunctionNet, 102
Mcs::Usb, 66	AreTransistorVoltagesSet, 104
ButterworthFilterBandPassNet, 99	CCMOSMea_FunctionNet, 104
ButterworthFilterBandPassNet, 99	ClearSTGOutput, 104
ButterworthFilterHighPassNet, 99	DetectChipType, 104
ButterworthFilterHighPassNet, 99	EnableChannelsInGroup, 104
ButterworthFilterLowPassNet, 100	GetADCInputOffset, 105
ButterworthFilterLowPassNet, 100	GetBath, 105
Butterworth mercowi assivet, 100	GetBathMode, 105
CalibrateThermocouple	GetEnabledChannelsInGroup, 105
CFluidControlDeviceNet, 144	GetGate, 105
CTcxDeviceNet, 589	GetGNDI, 105
Campden_Ci4600EphysVideoDataIntegrator	GetGroupADCBits, 105
Mcs::Usb, 77	GetGroupChannelBitmaskBySelect, 105, 106
CancelInternalCalibration	GetGroupChannelBitmaskHS1NCBathCurrent,
CTEERFunctionNet, 602	106
CancelPoolLoop	GetGroupChannelBitmaskHS1NCCol2Current,
CRoboDeviceNet, 474	106
CancelPoolLoopAndStopMovement	
CRoboDeviceNet, 474	GetGroupChannelBitmaskHS1NChipTemp, 106
CancelTableLoop	GetGroupChannelBitmaskHS1Sidebands, 106,
CRoboDacqNet, 459	107
CancelTableLoopAndStopTable	GetGroupChannelBitmaskHS1TriggerStatus, 107
CRoboDacqNet, 459	GetGroupChannelBitmaskIFDigChannels, 107
CapacityTest	GetGroupChannelBitmaskInterfaceADC, 107
CMultiBatteryChargerDeviceNet, 376	GetGroupChannelBitmaskPacketFrameContext,
CatchAmp	107, 108
Mcs::Usb, 75	GetGroupChannelBitmaskSTG1DACSignal, 108
	GetGroupDCOffset, 108
CatchAmpGetAdcMean	GetGroupID, 108
CMcsBus_SensorNet, 228	GetGroupNumberOfChannels, 108, 109
CatchAmpGetAdcValue	GetGroupResolutionPerDigit, 109
CMcsBus_SensorNet, 228	GetGroupSampleSize, 109
CatchAmpGetAdcValueH	GetGroupType, 109
CMcsBus_SensorNet, 228	GetGroupUnit, 109
CatchAmpGetAdcValueL	GetMaxNumOfColumns, 110
CMcsBus_SensorNet, 228	GetNeurochipMemoryData, 110
CatchAmpGetDacAmplitude	GetNeurochipMemorySize, 110
CMcsBus_SensorNet, 229	

GetNumberOfSupportedGroups, 110	WidthFrequency, 118
GetSourceBulk, 110	CDacCalibrationFunctionNet, 118
GetSourceDrain, 110	!CDacCalibrationFunctionNet, 119
GetSourceGate, 110	~CDacCalibrationFunctionNet, 119
GetStimulusSites, 110	BurnDacOffset, 119
GetVDD3I, 111	CDacCalibrationFunctionNet, 119
GetVDDI, 111	GetDacOffset, 119
IsChipPowered, 111	SetDacOffset, 120
IsGateFloating, 111	CDacqGroupChannelGenericSelectionNet, 120
PowerChip, 111	CDacqGroupChannelGenericSelectionNet, 120
SetADCInputOffset, 111	CDacqGroupChannelSelectionNet, 121
SetBath, 111	CDacqGroupChannelSelectionNet, 121
SetBathMode, 111	CDacqGroupChannelSelectionTemplateNet
SetGate, 111	CDacqGroupChannelSelectionTemplateNet< Dac-
SetGateFloating, 111	qGroupChannelEnumTemplateNet, Dac-
SetGateToVOP, 112	qGroupChannelEnumTemplate, CDevice-
SetNeurochipMemoryData, 112	GroupChannelInfoTemplateNet >, 122
SetSourceBulk, 112	CDacqGroupChannelSelectionTemplateNet< Dac-
SetSourceDrain, 112	qGroupChannelEnumTemplateNet, Dac-
SetSourceGate, 112	qGroupChannelEnumTemplate, CDevice-
SetStimulusSites, 112	GroupChannelInfoTemplateNet >, 121
UpdateTransistorVoltages, 112	CDacqGroupChannelSelectionTemplateNet, 122
VOPSTimerSetResetTimes, 112	EnableChannelsInGroup, 122
CCMOSMeaDeviceNet, 113	GetDeviceGroupChannelInfos, 122
~CCMOSMeaDeviceNet, 114	GetEnabledChannelsInGroup, 122, 123
CCMOSMeaDeviceNet, 114	GetGroupID, 123
CMosMea, 115	GetGroupNumberOfChannels, 123
GetAvailableBaseSamplerates, 114	GetGroupSampleSize, 123
GetBaseSamplerate, 114	GetGroupType, 123
GetChannelDatal16, 114	GetNumberOfSupportedGroups, 124
GetChannelDatal32, 114	CDeviceGroupChannelInfoGenericNet, 124
GetChannelDataUI16, 114	CDeviceGroupChannelInfoGenericNet, 124
GetChannelDataUl32, 115	CDeviceGroupChannelInfoMEA2100_256Net, 124
GetCMOSDataDictionary, 115	CDeviceGroupChannelInfoMEA2100_256Net, 125
GetMaxReadableColumns, 115	CDeviceGroupChannelInfoNet, 125
SetBaseSamplerate, 115	CDeviceGroupChannelInfoNet, 125
SetRegionOfInterests, 115	CDeviceGroupChannelInfoSCUNet, 126
Stimulus, 116	CDeviceGroupChannelInfoSCUNet, 126
UpdateChannelBlock, 115	CDeviceGroupChannelInfoTemplateNet
CCMOSMeaDeviceNet::CRegionOfInterestRect, 445	CDeviceGroupChannelInfoTemplateNet< Dacq-
CRegionOfInterestRect, 445	GroupChannelEnumTemplateNet >, 126
DeepCopy, 445	CDeviceGroupChannelInfoTemplateNet< Dacq-
m_Bottom, 446	GroupChannelEnumTemplateNet $>$, 126
m_Left, 446	CDeviceGroupChannelInfoTemplateNet, 126
m_Right, 446	GroupID, 127
m_Top, 446	GroupType, 127
CCreateFilterNet, 116	NumberOfChannels, 127
\sim CCreateFilterNet, 117	CDeviceGroupChannelInfoW2100Net, 127
CCreateFilterNet, 116, 117	CDeviceGroupChannelInfoW2100Net, 127
CutoffOrCenterFrequency, 118	CDigOutStimulatorFunctionNet, 128
FindFilter, 117	!CDigOutStimulatorFunctionNet, 129
GetBiQuad, 117	\sim CDigOutStimulatorFunctionNet, 129
GetBiQuads, 117	CDigOutStimulatorFunctionNet, 129
GetCpp, 117	ClearChannel, 129
NumCoefSets, 118	GetGlobalRepeat, 129
Order, 118	GetNumberOfChannels, 130
SampleRate, 118	GetStartTriggerSlope, 130
Scale, 118	GetStopTriggerSlope, 130

PrepareChannelData, 130	ToString, 142
SendChannelData, 131	CFirmwareDestinationNet
SetGlobalRepeat, 131	Mcs::Usb, 53
SetStartTriggerSlope, 131	CFluidControlDeviceNet, 143
SetStopTriggerSlope, 132	~CFluidControlDeviceNet, 144
CEncapsulatorDeviceNet, 132	CalibrateThermocouple, 144
CEncapsulatorDeviceNet, 132	CFluidControlDeviceNet, 144
GetRoboFluidDevice, 133	GetAdc, 144
CExternDTesterDeviceNet, 133	GetDigin, 145
!CExternDTesterDeviceNet, 134	GetDigout, 145
~CExternDTesterDeviceNet, 134	GetPWM, 145
CExternDTesterDeviceNet, 133	GetReferenceTemperature, 145
Read, 134	GetSingleValve, 146
Read2, 134	GetThermocoupleCalibration, 146
Write, 134	GetThermocoupleNanovoltPerKelvin, 146
Write2, 135	GetThermocoupleTemperature, 147
CFilterCoefficientsNet, 135	GetValve, 147
∼CFilterCoefficientsNet, 136	McsBus_VoltageMode, 149
A, 137	SetDigout, 147
B, 137	SetPWM, 147
CFilterCoefficientsNet, 135, 136	SetSingleValve, 148
GetUintA, 136	SetThermocoupleNanovoltPerKelvin, 148
GetUintB, 136	SetValve, 148
IsEqual, 136	CFYIDeviceNet, 149
CFilterCoefficientsNet::s_FilterAttributesNet, 709	CFYIDeviceNet, 149
CommaPositionA, 710	FYIProgram, 150
CommaPositionB, 710	FYITemp, 150
PostCommaA, 710	Sensor, 150
PostCommaB, 710	CGenericDevelopDeviceNet, 150
PreCommaA, 710	~CGenericDevelopDeviceNet, 157
PreCommaB, 710	CGenericDevelopDeviceNet, 157
s_FilterAttributesNet, 709	ClosePipe, 157
ToCpp, 709	FindEndpoints, 158
CFilterConfigurationNet, 137	GetBuffer, 158
CFilterConfigurationNet, 138	GetByteBuffer, 158
EraseFilterParameterPermanent, 138	GetIntBuffer, 158
GetDefaultFilterAttributes, 138	GetShortBuffer, 159
GetFilterAttributes, 138	GetUByteBuffer, 160
GetHighpassFilterEnable, 138	GetUIntBuffer, 160
ResetHighpassFilter, 138	GetUShortBuffer, 161
SetFilterParameter, 138, 139	OpenPipe, 162
SetFilterParameterPermanent, 139	ReadPipe, 162
SetHighpassFilterEnable, 139	ResetPipe, 163
CFilterConfigurationRegisterNet, 139	SetBuffer, 163
CFilterConfigurationRegisterNet, 140	SetByteBuffer, 163
EraseFilterParameterPermanent, 140	SetIntBuffer, 164
SetFilterParameter, 140	SetShortBuffer, 164
SetFilterParameterPermanent, 140, 141	SetUByteBuffer, 165
CFilterPropertyNet, 141	SetUIntBuffer, 165
~CFilterPropertyNet, 141	SetUShortBuffer, 166
CFilterPropertyNet, 141	SetValue, 167
CornerFrequency, 142	VendorInRequest, 167
CornerFrequencymHz, 142	VendorOutRequest, 167
FilterActive, 142	WritePipe, 167
FilterBand, 142	CGilsonDeviceNet, 168
FilterFamily, 142	~CGilsonDeviceNet, 169
FilterType, 142	CGilsonDeviceNet, 169
Order, 142	ConnectSlave, 169

CatlantAnguar 100	Channel Dlack Dead Frames 199
GetLastAnswer, 169	ChannelBlock_ReadFramesl32
m_pGilsonDevice, 169	CMcsUsbDacqNet, 266
SendBuffered, 169	ChannelBlock_ReadFramesUI16
SendImmediate, 169	CMcsUsbDacqNet, 267
SendImmediateGetResponse, 169	ChannelBlock_ReadFramesUl32
CGrapheneASICDeviceNet, 170	CMcsUsbDacqNet, 268, 269
\sim CGrapheneASICDeviceNet, 170	channeldata_current
CGrapheneASICDeviceNet, 170	Mcs::Usb, 85
GetAvailableBaseSamplerates, 170	channeldata_current_always_boost
SetBaseSamplerate, 171	Mcs::Usb, 86
SetRegionOfInterests, 171	channeldata_current_always_boost_own_sync
CGrapheneFunctionNet, 171	Mcs::Usb, 86
!CGrapheneFunctionNet, 173	channeldata_current_own_boost_gnd_sync
~CGrapheneFunctionNet, 173	Mcs::Usb, 86
CGrapheneFunctionNet, 173	channeldata_current_own_sync
GetADCGain, 173, 174	Mcs::Usb, 86
GetADCOffset, 174	channeldata_positive_current
GetCur2VolResistance, 175	Mcs::Usb, 86
GetDACOffset, 175	channeldata_positive_current_own_boost_gnd_sync
GetVdsVgs, 176	Mcs::Usb, 86
GetVdVs, 176, 177	channeldata_positive_current_own_sync
GetVdVsDAC, 177	Mcs::Usb, 86
GetVoltageRange, 177, 178	channeldata_positive_voltage
GetVoltageReached, 178	Mcs::Usb, 86
GetVoltageResolution, 178, 179	
	channeldata_voltage
SetADCGain, 179	Mcs::Usb, 85
SetADCOffset, 180	ChannelDataEvent
SetDACOffset, 180	CMcsUsbDacqNet, 297
SetVdsVgs, 181	ChannelPIC
SetVdVs, 181	Mcs::Usb, 55
SetVdVsDAC, 182	ChannelReset
SetVoltageRange, 182, 183	CMultiBatteryChargerDeviceNet, 376
SetVoltageResolution, 183	ChannelTest
ChangeSerialNumber	Mcs::Usb, 78
CMcsUsbFactoryNet, 301	Charge
channel	BatteryState, 96
CMeaAudioFunctionNet::s_setaudionet, 710	ChargeRegionString
CW2100_FunctionNet::AudioChannelsNet, 96	BatteryState, 96
ChannelBlock_AvailFrames	ChargeString
CMcsUsbDacqNet, 256	BatteryState, 96
ChannelBlock_ReadAsFrameArrayI16	ChecksumAndPacketCounter
CMcsUsbDacqNet, 256, 257	Mcs::Usb, 60
ChannelBlock_ReadAsFrameArrayl32	CHiClampDeviceNet, 183
CMcsUsbDacqNet, 258, 259	CHiClampDeviceNet, 184
ChannelBlock_ReadAsFrameArrayUI16	RoboDacq, 184
CMcsUsbDacqNet, 259, 260	CHLADacqNet, 184
ChannelBlock ReadAsFrameArrayUI32	CHLADacqNet, 185
CMcsUsbDacqNet, 261	CHLADeviceNet, 185
ChannelBlock ReadFramesDictI16	CHLADeviceNet, 185
CMcsUsbDacqNet, 262	HLADacq, 186
ChannelBlock_ReadFramesDictl32	SerialPort, 186
CMcsUsbDacqNet, 263	CHWInfo
ChannelBlock ReadFramesDictUI16	CMcsUsbDacqNet::CHWInfo, 187
CMcsUsbDacqNet, 263	Ci4600Intan
ChannelBlock ReadFramesDictUI32	Mcs::Usb, 52
CMcsUsbDacqNet, 264	CIntanMea_FunctionNet, 189
ChannelBlock ReadFramesI16	AmplifierSettle, 190
CMcsUsbDacqNet, 264, 265	BeginImpedanceCheck, 190
01v103030Da0q1vet, 204, 203	Degininpedanceoneck, 190

CIntanMea_FunctionNet, 189	ClearUserDefinedNameCache
GetDSPHighPassByIndex, 190	CW2100 FunctionNet, 616
	ClearValveTable
GetImpedanceResult, 190	
GetIntanRegister, 190	CWarnerValveControllerDeviceNet, 651
GetLowerFrequencyByIndex, 190	CLIH3DeviceNet, 195
GetUpperFrequencyByIndex, 190	!CLIH3DeviceNet, 197
SetBandwidthByIndex, 190	∼CLIH3DeviceNet, 197
SetDiagnosticMode, 190	CLIH3DeviceNet, 197
SetDSPHighPassByIndex, 191	DummyCommand, 197
SetIntanRegister, 191	EnableUserTrigger, 197
CInterfaceboard2FunctionNet, 191	ErasePermanentAdcOffset, 198
!CInterfaceboard2FunctionNet, 192	ErasePermanentDacOffset, 198
\sim CInterfaceboard2FunctionNet, 192	GetAdcOffset, 198
CInterfaceboard2FunctionNet, 192	GetAudioOutDacParameter, 198
GetloVoltage, 192	GetDacIdleValue, 199
SetloVoltage, 192	GetDacOffset, 199
CInterfaceboardFunctionNet, 193	GetDacqRunStatus, 199
!CInterfaceboardFunctionNet, 194	GetDacUseIdleValue, 199
~CInterfaceboardFunctionNet, 194	GetDigInState, 200
CInterfaceboardFunctionNet, 193	GetEEpromPage, 200
GetCardinalDacqSamplerate, 194	GetSampleInterval, 200
GetCardinalStgOutputrate, 194	IsUserTriggerEnabled, 200
SetCardinalDacqSamplerate, 194	ReadClipping, 201
SetCardinalStgOutputrate, 194	ReadUARTData, 201
ClampAmpRestart	SendCommand, 201
CRoboDacqNet, 459	SetAdcOffset, 201
ClampModeCurrent	SetAdcOffsetPermanent, 202
Mcs::Usb, 87	SetAudioOutDacParameter, 202
ClampModeInternalCalibration	SetDacIdleValue, 202
Mcs::Usb, 87	SetDacOffset, 202
ClampModeOpen	SetDacOffsetPermanent, 203
Mcs::Usb, 87	SetDacUseIdleValue, 203
ClampModeVoltage	SetDigOutState, 203
Mcs::Usb, 87	SetEEpromPage, 203
ClearBuffers	SetSampleInterval, 204
CMcsUsbDacqNet, 269	StimulusFunction, 204
ClearChannel	WriteUARTData, 204
CDigOutStimulatorFunctionNet, 129	Close
ClearChannel PrepareAndSendData	Mcs::Usb, 76
CStg200xDownloadNet, 568	CloseAllValves
CStimulusFunctionNet, 575	CRoboFluidDeviceNet, 488
ClearChannelData	ClosePipe
	•
CStg200xDownloadBasicNet, 561	CGenericDevelopDeviceNet, 157 ClosePlateClamp
CStimulusFunctionNet, 575	•
CW2100_StimulatorFunctionNet, 623	CMultiwellDeviceNet, 386
ClearMultiplexedData	CMcsBus_AxisParametersNet, 204
CStimulusFunctionNet, 575	~CMcsBus_AxisParametersNet, 205
ClearSTGOutput	CMcsBus_AxisParametersNet, 205
CCMOSMea_FunctionNet, 104	GetAxisParametersSignedEeprom, 205
ClearStimulusParametersCache	GetAxisParametersUnsignedEeprom, 205
CW2100_FunctionNet, 616	SetAxisParametersEeprom, 205, 206
ClearSyncData	CMcsBus_ExtensionNet, 206
CStg200xDownloadBasicNet, 561	~CMcsBus_ExtensionNet, 206
CStimulusFunctionNet, 575	CMcsBus_ExtensionNet, 206
ClearTable	GetLEDSwitch, 207
CRetinaLedDeviceNet, 447	SetLEDSwitch, 207
ClearTableName	CMcsBus_FYIExtensionNet, 207
CWarnerValveControllerDeviceNet, 651	~CMcsBus FYIExtensionNet. 207

CMcsBus_FYIExtensionNet, 207	SetMCAccelerationShortCommand, 219
GetDIO, 208	SetMCAxisRevisionEeprom, 219
GetSingleHeater, 208	SetMCBreakCurrent, 219
GetValves, 208	SetMCBreakCurrentEeprom, 219
SetDIO, 208	SetMCConfig, 219
SetSingleHeater, 208	SetMCConfigEeprom, 220
SetValves, 208	SetMCCurrent, 220
CMcsBus_MotorControlNet, 209	SetMCCurrentEeprom, 220
~CMcsBus_MotorControlNet, 212	SetMCCurrentMode, 220
CMcsBus MotorControlNet, 212	SetMCCurrentModeEeprom, 220
GetMCAcceleration, 212	SetMCCurrentModeShortCommand, 220
GetMCAccelerationEeprom, 212	SetMCCurrentPosition, 221
GetMCAccelerationShortCommand, 212	SetMCCurrentShortCommand, 221
GetMCAxisRevisionEeprom, 212	SetMCMaxAcceleration, 221
GetMCBreakCurrent, 212	SetMCMaxAccelerationEeprom, 221
GetMCBreakCurrentEeprom, 213	SetMCMaxCurrent, 221
GetMCConfig, 213	SetMCMaxCurrentEeprom, 221
GetMCConfigEeprom, 213	SetMCMaxSpeed, 222
GetMCCurrent, 213	SetMCMaxSpeedEeprom, 222
GetMCCurrentEeprom, 213	SetMCMaxTravel, 222
GetMCCurrentMode, 213	SetMCMaxTravelEeprom, 222
GetMCCurrentModeEeprom, 213	SetMCMaxTravelShortCommand, 222
GetMCCurrentModeShortCommand, 214	SetMCNewPosition, 222
GetMCCurrentPosition, 214	SetMCOutputOnOff, 223
GetMCCurrentShortCommand, 214	SetMCReference, 223
GetMCCurrentSpeed, 214	SetMCReferenceCurrent, 223
GetMCMaxAcceleration, 214	SetMCReferenceCurrentEeprom, 223
GetMCMaxAccelerationEeprom, 214	SetMCRegulatorGain, 223
GetMCMaxCurrent, 214	SetMCRegulatorGainEeprom, 223
GetMCMaxCurrentEeprom, 215	SetMCRotation, 224
GetMCMaxSpeed, 215	SetMCScalingFactor, 224
GetMCMaxSpeedEeprom, 215	SetMCScalingFactorEeprom, 224
GetMCMaxTravel, 215	SetMCSpeed, 224
GetMCMaxTravelEeprom, 215	SetMCSpeedEeprom, 224
GetMCMaxTravelShortCommand, 215	SetMCSpeedShortCommand, 224
GetMCMovement, 215	SetMCSpeedUnitEeprom, 225
GetMCNewPosition, 216	SetMCStandbyCurrent, 225
GetMCOutputOnOff, 216	SetMCStandbyCurrentEeprom, 225
GetMCPhase, 216	SetMCStandbyTime, 225
GetMCPhaseOffset, 216	SetMCStandbyTimeEeprom, 225
GetMCReference, 216	SetSubChannel, 225
GetMCReferenceCurrent, 216	StartMCMovement, 226
GetMCReferenceCurrentEeprom, 216	StopMCMovement, 226
GetMCRegulatorGain, 217	CMcsBus SensorNet, 226
GetMCRegulatorGainEeprom, 217	~CMcsBus_SensorNet, 228
GetMCScalingFactor, 217	CatchAmpGetAdcMean, 228
GetMCScalingFactorEeprom, 217	CatchAmpGetAdcValue, 228
GetMCSpeed, 217	CatchAmpGetAdcValueH, 228
GetMCSpeedEeprom, 217	CatchAmpGetAdcValueL, 228
GetMCSpeedShortCommand, 217	CatchAmpGetDacAmplitude, 229
GetMCSpeedUnitEeprom, 218	CatchAmpGetDacEnable, 229
GetMCStandbyCurrent, 218	CatchAmpGetDacOffset, 229
GetMCStandbyCurrentEeprom, 218	CatchAmpGetPwmEnable, 229
GetMCStandbyTime, 218	CatchAmpSetDacAmplitude, 229
GetMCStandbyTimeEeprom, 218	CatchAmpSetDacEnable, 229
GetSubChannel, 218	CatchAmpSetDacOffset, 229
SetMCAcceleration, 218	CatchAmpSetPwmEnable, 229
SetMCAccelerationEeprom, 219	CMcsBus_SensorNet, 228

Get2AnalogInput, 230	\sim CMcsBus_VoltageModeNet, 239
Get2DigitalInput, 230	CMcsBus VoltageModeNet, 239
Get4ADC, 230	GetVMMaxNegativeCurrent, 240
Get4ADCAverage, 230	GetVMMaxNegativeCurrentEeprom, 240
Get4ADCCatchampAverageShift, 230	GetVMMaxNegativeVoltage, 240
Get4ADCMode, 230	GetVMMaxNegativeVoltageEeprom, 240
Get4DAC, 230	GetVMMaxPositiveCurrent, 240
GetADCs, 230	GetVMMaxPositiveCurrentEeprom, 240
GetADCsLoop, 231	GetVMMaxPositiveVoltage, 240
GetBubbleStatus, 231	GetVMMaxPositiveVoltageEeprom, 241
GetDACs, 231	GetVMOutputOnOff, 241
GetDAO3, 201 GetDetectionThreshold, 231	GetVMVoltage, 241
GetDetectorValue, 231	SetVMMaxNegativeCurrent, 241
GetLatency, 231	SetVMMaxNegativeCurrentEeprom, 241
GetLatencyCounter, 231	SetVMMaxNegativeVoltage, 241
GetMinimalThreshold, 231	SetVMMaxNegativeVoltageEeprom, 241
GetMovePump, 232	SetVMMaxPositiveCurrent, 242
GetPiezoState, 232	SetVMMaxPositiveCurrentEeprom, 242
GetPressure, 232	SetVMMaxPositiveVoltage, 242
GetPressureOffset, 232 GetRegulationTimeouts, 232	SetVMMaxPositiveVoltageEeprom, 242 SetVMOutputOnOff, 242
GetRegulatorFactor, 233	SetVMoutputOrion, 242 SetVMVoltage, 242
GetRegulatorOnOff, 233	CMcsBusNet, 243
GetRegulatorStatus, 233	
GetRotatePump, 233	~CMcsBusNet, 243
•	CMcsBusNet; 243
GetSamplePeriode, 233	CMcsBusNet::GetMode, 244
GetSuneState 223	CMcsBusNet::GetModeEeprom, 244
GetSyncState, 233	CMcsBusNet::SetMode, 244
Set4ADCCatchampAverageShift, 234	CMcsBusNet::SetModeEeprom, 244
Set4ADCMode, 234	GetBusAddress, 244
Set4DAC, 234	GetBusAddressEeprom, 244
SetDACs, 234	Get IV/Povision Forgross 045
SetDetectionThreshold, 234	GetHWRevisionEeprom, 245
SetLatency, 234	SetBusAddress, 245
SetMinimalThreshold, 234	SetBusAddressEeprom, 245
SetMovePump, 235	SetCommand, 246
SetPiezoState, 235	SetHWRevisionEeprom, 246
SetPressureOffset, 235	CMcsBusNet::GetMode
SetRegulationTimeouts, 235	CMcsBusNet, 244
SetRegulatorFactor, 235	CMcsBusNet::GetModeEeprom
SetRegulatorOnOff, 235	CMcsBusNet, 244
SetRotatePump, 235	CMcsBusNet::SetMode
SetSamplePeriode, 236	CMcsBusNet, 244
SetSollPressure, 236	CMcsBusNet::SetModeEeprom
StartSync, 236	CMcsBusNet, 244
TactSwitchGetState, 236	CMcsUsbDacqNet, 247
TactSwitchSetDisplay, 236	~CMcsUsbDacqNet, 253
CMcsBus_TempSensorNet, 236	AddSelectedChannelsQueue, 253–256
~CMcsBus_TempSensorNet, 237	ChannelBlock_AvailFrames, 256
CMcsBus_TempSensorNet, 237	ChannelBlock_ReadAsFrameArrayl16, 256, 257
GetNanoVoltsPerKelvin, 237	ChannelBlock_ReadAsFrameArrayl32, 258, 259
GetTemperatur, 237	ChannelBlock_ReadAsFrameArrayUI16, 259, 260
GetThermoOffset, 237	ChannelBlock_ReadAsFrameArrayUl32, 261
GetThermoTemp, 238	ChannelBlock_ReadFramesDictl16, 262
GetThermoVoltage, 238	ChannelBlock_ReadFramesDictl32, 263
SetNanoVoltsPerKelvin, 238	ChannelBlock_ReadFramesDictUI16, 263
SetThermoOffset, 238	ChannelBlock_ReadFramesDictUl32, 264
CMcsBus_VoltageModeNet, 238	ChannelBlock_ReadFramesI16, 264, 265

ChannelBlock_ReadFramesl32, 266	SetVoltageRangeInMicroVolt, 292
ChannelBlock_ReadFramesUI16, 267	StartDacq, 292, 293
ChannelBlock_ReadFramesUl32, 268, 269	StartLoop, 294, 295
ChannelDataEvent, 297	StopDacq, 296
ClearBuffers, 269	StopLoop, 296
CMcsUsbDacqNet, 253	CMcsUsbDacqNet::CHWInfo, 186
CMcsUsbDacqNet::GetFilterProperties, 269	CHWInfo, 187
Error_Callback_Aquisition_Stopped, 296	GetAvailableSampleRates, 187
Error_Callback_Data_lost, 296	GetAvailableVoltageRangesInMicroVolt, 187
Error_Callback_Frames_Lost, 296	GetAvailableVoltageRangesInMicroVoltAnd-
Error_Callback_Packet_Error, 296	StringsInMilliVolt, 187
Error_Callback_Queue_Full, 297	GetNumberOfHWADCChannels, 188
Error_Callback_RingQueue_Full, 297	GetNumberOfHWDigitalChannels, 188
ErrorEvent, 297	IsDigitalChannelDedicated, 188
GetAdapterType, 269	CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet
GetAdcDataFormat, 270	613
GetAdcZero, 270	CVoltageRangeInfoNet, 614
GetAnalogValueUnit, 270	VoltageRangeDisplayStringMilliVolt, 614
GetChannelDataFillSize, 270	VoltageRangeInMicroVolt, 614
GetChannelLayout, 270	CMcsUsbDacqNet::GetFilterProperties
GetChannelsInBlock, 270	CMcsUsbDacqNet, 269
GetDataFormat, 271	CMcsUsbDeviceStatePushFunctionNet, 297
GetDataMode, 271	CMcsUsbDeviceStatePushFunctionNet, 298
GetDigitalSource, 271–273	McsUsbDeviceStateEvent, 298
GetErrorMessage, 274	TriggerStatus, 298
GetFilterProperty, 274	CMcsUsbDeviceStatePushNet, 298
GetGroupChannelDataI16, 274	CMcsUsbDeviceStatePushNet, 299
GetGroupChannelDatal32, 275	McsUsbDeviceStateEvent, 299
GetGroupChannelDataUI16, 275	TriggerStatus, 299
GetGroupChannelDataUI32, 276	CMcsUsbFactoryNet, 299
GetHardwareMaxRange, 276	~CMcsUsbFactoryNet, 301
GetHardwareMinRange, 276	ChangeSerialNumber, 301
GetMaxSamplingFrequency, 277	CMcsUsbFactoryNet, 301
GetMeaLayout, 277	Coldstart, 301
GetMinSamplingFrequencyStepsize, 277	CompareFirmware, 302
GetNumberOfDataBits, 277	DownloadFirmware, 302
GetPoti, 277	FindFirmwareVersionMagicInBuffer, 302
GetResolutionPerDigit, 278	FX3MCSDataAddress, 307
GetSamplerate, 278	FX3MCSDataDeviceIdOffset, 307
GetVoltageRangeIndex, 278	FX3MCSDataIFB1ImageOffset, 307
GetVoltageRangeInMicroVolt, 278	FX3MCSDataIFB2ImageOffset, 307
GetVoltageRangeInMilliVolt, 278	FX3MCSDataVersionOffset, 307
HWInfo, 278	GetChecksumFromFX3Image, 302
Samplerate, 297	GetDestination, 302
SendStartDacq, 279	GetDestinationDisplayLabel, 302
SendStartStgAndDacq, 279	GetDestinationName, 302, 303
SendStopDacq, 279	GetDestinationSerialNumber, 303
SendStopStgAndDacq, 280	GetDestinationTargetAddress, 303
SendStopStgAndDacqWithOptions, 280	GetDestinationTargetSelectionInfo, 303
SetDataMode, 280	GetFirmwareVersionFromFile, 303
SetDigitalSource, 281–283	GetFirmwareVersionFromHexFile, 303
SetPoti, 283	GetNumDestinations, 304
SetSamplerate, 283	GetUSBDeviceIDFromFX3Image, 304
SetSelectedChannels, 285–287	GetUsercodeFromBitFile, 304
SetSelectedChannelsQueue, 287–289	GetUsercodeFromFlash, 304
SetSelectedData, 289–291	GetXilinxFlashOffset, 304
SetupGroupDacqQueue, 291	GetXilinxFlashReadCommand, 304
SetVoltageRangeByIndex, 291	LoadUserFirmware, 304, 305

5 15 15 5 5 1	
ReadBlockFromFlash, 305	GetDeviceRootHubVendorEnum, 325
ReadBlockFromIFBGlobalEEprom, 305	GetDeviceRootHubVendorID, 325
ReadBlockFromNVMEM, 305	GetDeviceRootHubVendorName, 325
SetDestinationSerialNumber, 305	GetDeviceSpeed, 326
UpdateFirmware, 305–307	GetErrorText, 326
CMcsUsbFactoryNet::TargetSelectionInfoNet, 714	GetFirmwareVersion, 326
Address, 714	GetHardwareRevision, 326
Value, 714	GetHeadstageActive, 327
CMcsUsbFunctionNet, 308	GetHeadstageID, 327
!CMcsUsbFunctionNet, 309	GetHeadstagePresent, 327
~CMcsUsbFunctionNet, 309	GetLastUSBError, 327
CMcsUsbFunctionNet, 308, 309	GetMea21UsbPort, 328
m_pMcsUsb, 309	GetNumConfigurations, 328
m_pMcsUsbFunction, 309	GetSerialNumber, 328
ThrowCUsbExceptionNetOnError, 309	GetSoftwareKey, 328
CMcsUsbFunctionPointerContainer, 309	GetSoftwareKeyString, 328
CMcsUsbListEntryNet, 309	GetStatus, 328
~CMcsUsbListEntryNet, 310	GetStatusOfLastCommand, 328
Deviceld, 314	GetUsbListEntry, 329
DeviceName, 314	GetVersion, 329
Equals, 311	HasSoftwareKey, 329
GetEntry, 312	• •
• •	IsConnected, 329
GetEntryCount, 313	IsDeviceHighSpeed, 329
HwVersion, 314	IsDeviceHighSpeedCapable, 329
Manufacturer, 314	IsExceptionsEnabled, 330
Product, 314	MultibootGetCypressImageId, 330
SerialNumber, 314	MultibootGetImageId, 330
SetStringFormat, 313	MultibootGetSelectedImage, 330
ToString, 314	MultibootSelectImage, 330
CMcsUsbListNet, 315	ReadEepromRegisterPreconfig, 330, 331
!CMcsUsbListNet, 316	ReadRegister, 331
\sim CMcsUsbListNet, 316	ReadRegister32, 331
CMcsUsbListNet, 315	ReadRegisterTimeSlot, 331
Count, 317	RemoveSoftwareKey, 331
DeviceArrival, 317	RescanHeadstage, 331
DeviceRemoval, 317	SerialNumber, 338
GetNumberOfDevices, 316	SetConfiguration, 332
GetUsbListEntries, 316	SetSoftwareKey, 332
GetUsbListEntry, 316	Status_AlreadyConfigured, 334
IsDeviceTypeOf, 317	Status BadStartFrame, 334
SetStringFormat, 317	Status_Btstuff, 334
CMcsUsbNet, 318	Status BufferOverrun, 334
!CMcsUsbNet, 322	Status BufferUnderrun, 334
∼CMcsUsbNet, 322	Status Canceled, 334
AddSoftwareKey, 322	Status_Canceling, 335
AssociateToThis, 322	Status ConnectedPipes, 335
CMcsUsbNet, 322	Status_ControlNotOwned, 335
Connect, 322, 323	Status_Crc, 335
CyclePort, 324	Status DataOverrun, 335
Disconnect, 324	Status_DataToggleMismatch, 335
EmptyKey, 324	Status_DataToggleMismatch, 333
EnableExceptions, 324	Status_DataOnderrun, 335 Status_DeviceLocked, 335
·	
EraseEepromRegisterPreconfig, 324	Status_DeviceNotFound, 335
GetConfiguration, 325	Status_DeviceRemoved, 335
GetDeviceCannotStallOutRequests, 325	Status_DevNotResponding, 335
GetDeviceCapableSpeed, 325	Status_EndpointHalted, 336
GetDeviceEnum, 325	Status_ErrorBusy, 336
GetDeviceId, 325	Status_ErrorShortTransfer, 336

Status_Fifo, 336	\sim CMeaCleanDeviceNet, 345
Status_FrameControlOwned, 336	CMeaCleanDeviceNet, 345
Status_InternalHcError, 336	GetCycle, 345
Status_InvalidDeviceHandle, 336	GetCycles, 345
Status_InvalidHandle, 336	GetMaxVoltage, 345
Status InvalidParameter, 336	GetMinVoltage, 346
Status_InvalidPipeHandle, 336	GetOutputVoltage, 346
Status InvalidUrbFunction, 337	GetSlope, 346
Status_loPending, 337	IsRunning, 346
Status IoTimeout, 337	SetCycles, 346
Status IsochRequestFailed, 337	SetMaxVoltage, 347
Status LastUsbErrorMismatch, 337	SetMinVoltage, 347
Status_NoBandwidth, 337	_
	SetSlope, 347
Status_NoMemory, 337	Start, 347
Status_NoSuchDevice, 337	Stop, 347
Status_NotAccessed, 337	CMeaCoatDeviceNet, 348
Status_NotSupported, 337	!CMeaCoatDeviceNet, 349
Status_PidCheckFailure, 337	~CMeaCoatDeviceNet, 349
Status_PipeNotLinked, 338	CMeaCoatDeviceNet, 349
Status_RequestFailed, 338	GetCurrentCycle, 349
Status_RequestMutexFailed, 338	GetCycles, 349
Status_RequestMutexTimeout, 338	GetDuration, 350
Status_Stall, 338	GetMaxCurrent, 350
Status_Unconfigured, 338	GetOffsetCurrent, 350
Status_UnexpectedPid, 338	GetOutputCurrent, 350
ThrowCUsbExceptionNetOnError, 332	GetPauseDuration, 350
TxnGetSerialNumber, 332	GetSlope, 351
TxnSetSerialNumber, 332	GetTimeInPause, 351
TxnTestMemoryReadAndCheck, 332	GetTimeInPlateau, 351
TxnTestMemoryWrite, 332	IsRunning, 351
ValidKey, 332	SetCycles, 351
WPAError_ScanningIsPending, 338	SetDuration, 352
WriteEepromRegisterPreconfig, 333	SetMaxCurrent, 352
WriteRegister, 333	SetOffsetCurrent, 352
WriteRegister32, 333	SetPauseDuration, 352
WriteRegisterArray, 333	SetSlope, 353
WriteRegisterTimeSlot, 334	Start, 353
WriteRegisterValue, 334	Stop, 353
CMcsUsbPointerContainer, 339	CMeaDeviceNet, 353
CMEA2100_256DacqGroupChannelSelectionNet, 339	~CMeaDeviceNet, 355
CMEA2100 256DacqGroupChannelSelectionNet,	AnalogGain, 360
339	CMeaDeviceNet, 355
CMEA2100x256FunctionNet, 339	EnableChecksum, 355
!CMEA2100x256FunctionNet, 340	EnableDigitalIn, 356
~CMEA2100x256FunctionNet, 340	-
CMEA2100x256FunctionNet, 340	EnableTimestamp, 357
•	Gain, 360
GetLayoutConfiguration, 340	GetAnalogGain, 357
SetLayoutConfiguration, 340	GetEnumerationSpeed, 357
CMeaAudioFunctionNet, 341	GetGain, 357
CMeaAudioFunctionNet, 341	MeaAudioFunctionNet, 361
GetAudioChannels, 342	MeaDigitalDataFunctionNet, 361
GetNumberOfAudioChannels, 342	MeaFeedbackFunctionNet, 361
SetAudioChannels, 343	MeFunctionNet, 361
CMeaAudioFunctionNet::s_setaudionet, 710	SetDigitalOut, 357
amplification, 710	SetNumberOfAnalogChannels, 358
channel, 710	SetNumberOfChannels, 359
CMeaCleanDeviceNet, 344	SetTriggerMaskValue, 359
!CMeaCleanDeviceNet, 345	SetTriggerPeriod, 360

W2100_FunctionNet, 361	!CMeFunctionNet, 374
WClassicFunctionNet, 361	∼CMeFunctionNet, 374
CMeaDigitalDataFunctionNet, 361	CMeFunctionNet, 374
CMeaDigitalDataFunctionNet, 362	SetTransformer, 374
GetDigitalData, 362	CMosMea
SetDigitalData, 362, 363	CCMOSMeaDeviceNet, 115
CMeaFeedbackFunctionNet, 363	CmosMea
CMeaFeedbackFunctionNet, 364	Mcs::Usb, 68
FeedbackGetSampleTimerCount, 364	CMOSMEA5000DigitalSourceEnumNet
FeedbackGetTriggerTotzeitFactor, 364	Mcs::Usb, 57
FeedbackSetAnalogSource, 365	CMOSMeaBathModeEnumNet
FeedbackSetChannelFilter, 365	Mcs::Usb, 57
FeedbackSetCreateFilter, 365	CmosMeaHeadstage
FeedbackSetDigitalMapping, 365	Mcs::Usb, 67
FeedbackSetFeedback, 365	CMOSMeaHeadstage1NCBathCurrentEnumNet
FeedbackSetFilterOff, 365	Mcs::Usb, 57
FeedbackSetFilterParameter, 365	CMOSMeaHeadstage1NCCol2CurrentEnumNet
FeedbackSetFilterParameter32, 365	Mcs::Usb, 58
FeedbackSetGlobalChannelFilter, 366	CMOSMeaHeadstage1NChipTempEnumNet
FeedbackSetIIRFilterParameter, 366	Mcs::Usb, 58
FeedbackSetLogic, 366	CMOSMeaHS1SidebandEnumNet
FeedbackSetMkFilter, 366	Mcs::Usb, 58
FeedbackSetNumberOfLogics, 366	CMOSMeaHS1TriggerStatusEnumNet
FeedbackSetNumberOfRateCounter, 366	Mcs::Usb, 58
FeedbackSetNumberOfRateDetectors, 366	CmosmealFB2
FeedbackSetNumberOfSpikeDetectors, 367	Mcs::Usb, 67
FeedbackSetNumberOfTriggers, 367	CMOSMealFDigChannelEnumNet
FeedbackSetRateCounter, 367	Mcs::Usb, 59
FeedbackSetRateDetector, 367	CMOSMeaInterfaceADCEnumNet
FeedbackSetSpikeDetectorThreshold, 367	Mcs::Usb, 59
FeedbackSetTrigger, 367	CmosMeaInterfaceboard
FeedbackSetTriggerTotzeitFactor, 367	Mcs::Usb, 67
CMealmpedanceDeviceNet, 368	CMOSMeaPacketFrameContextGroupEnumNet
~CMeaImpedanceDeviceNet, 368	Mcs::Usb, 59
CMealmpedanceDeviceNet, 368	CMOSMeaSTG1DACSignalEnumNet
GetAdapterCode, 368	Mcs::Usb, 60
GetArraySize, 368	CMOSMeaValueUnitEnumNet
GetImpedanceTestFrequency, 369	Mcs::Usb, 60
GetReady, 369	CMultiBatteryChargerDeviceNet, 374
GetResult, 369	!CMultiBatteryChargerDeviceNet, 376
SetImpedanceTestFrequency, 369	~CMultiBatteryChargerDeviceNet, 376
StartMeasurement, 369	· · · · · · · · · · · · · · · · · · ·
	CapacityTest, 376
CMeasure Table DeviceNet, 369	ChannelReset, 376
CMeasureTableDeviceNet, 370	CMultiBatteryChargerDeviceNet, 376
Sensor, 370	GetBatteryVoltage, 377
CMeaSwitchDeviceNet, 370	GetChannels, 377
~CMeaSwitchDeviceNet, 371	GetChannelState, 377
CMeaSwitchDeviceNet, 371	GetChargeCapacity, 377
GetNumber, 371	GetChargeCurrent, 378
GetPattern, 371	GetChargingMode, 378
GetPatternBool, 371	GetChargingPCoefficient, 378
SetPattern, 371	GetDischargeCapacity, 378
SetPatternBool, 372	GetDischargeCurrent, 379
CMeaUSBDeviceNet, 372	GetDischargeCurrentSetPoint, 379
~CMeaUSBDeviceNet, 373	GetFinalDischargeVoltage, 379
CMeaUSBDeviceNet, 372, 373	GetRatedCapacity, 380
SetVoltageRangeForIndex, 373	SetChargingMode, 380
CMeFunctionNet, 373	SetChargingPCoefficient, 380

SetDischargeCurrentSetPoint, 380	COctoPotDeviceNet, 398
SetFinalDischargeVoltage, 381	EnableChecksum, 398
SetRatedCapacity, 381	EnableDigitaIIn, 398
SetRatedCapacityVolatile, 381	EnableTimestamp, 398
CMultiwellCallbackFunctionNet, 382	GetAdcOffset, 398
!CMultiwellCallbackFunctionNet, 383	GetDacOffset, 398
~CMultiwellCallbackFunctionNet, 383	PatternListStart, 399
CMultiwellCallbackFunctionNet, 382	RampStart, 399
GetPlateClampStateByHeadstage, 383	ResetAdcOffset, 399
GetPlateClampStateByHeadstageEvent, 383	ResetDacOffset, 399
OnGetPlateClampStateByHeadstage, 383	SetAdcOffset, 399
CMultiwellDeviceNet, 384	SetAmplificationSwitch, 399
!CMultiwellDeviceNet, 385	SetBathclamp, 399
~CMultiwellDeviceNet, 385	SetChannelSwitch, 399
ClosePlateClamp, 386	SetDacAutoControl, 399
CMultiwellDeviceNet, 385	SetDacOffset, 400
GetPlateClampLockState, 386	SetDacValue, 400
GetPlateClampState, 386	SetNumberOfChannels, 400
GetPlateMux, 386, 387	SetOutputRate, 400
GetPlateType, 387	SetPatternListEntry, 400
GetPowerMuxPlate, 387	SetPidParameter, 400
GetTouchPadEnable, 388	SetRampParameter, 400
GetVolatileClampOffset, 388	SetSineParameter, 400
IsPlateTypeValid, 388	SineStart, 401
LockPlateClamp, 389	
•	COkuvisionStimulatorDeviceNet, 401
OpenPlateClamp, 389	~COkuvisionStimulatorDeviceNet, 402
SetPlateMux, 389	COkuvisionStimulatorDeviceNet, 402
SetPlateType, 390	GetCheckVoltage, 402
SetPowerMuxPlate, 390	GetCurrentFactor, 402
SetTouchPadEnable, 390	GetDACOffset, 402
SetVolatileClampOffset, 391	GetMaxPower, 402
StopPlateClamp, 391	GetMaxVoltage, 402
UnlockPlateClamp, 391	GetPulseform, 402
CMultiwellOptoStimFunctionNet, 391	GetRTC, 403
!CMultiwellOptoStimFunctionNet, 392	GetStimulatorStatus, 403
\sim CMultiwellOptoStimFunctionNet, 392	GetVoltage, 403
CMultiwellOptoStimFunctionNet, 392	SetCheckVoltage, 403
GetAbsMaxCurrentInMicroAmp, 392	SetCurrentFactor, 403
GetColorRgb, 393	SetDACOffset, 403
GetColorStr, 393	SetMaxPower, 404
GetMaxDurationHighCurrentInMicroSec, 393	SetMaxVoltage, 404
GetMaxDutyCycleHighCurrent, 393	SetPulseform, 404
GetPermanentCurrentInMicroAmp, 394	SetRTC, 404
GetWaveLengthInNanometer, 394	Coldstart
SetAbsMaxCurrentInMicroAmp, 394	CMcsUsbFactoryNet, 301
SetColorRgb, 395	CommaPositionA
SetColorStr, 395	CFilterCoefficientsNet::s FilterAttributesNet, 710
SetMaxDurationHighCurrentInMicroSec, 395	Mcs::Usb, 65
SetMaxDutyCycleHighCurrent, 395	CommaPositionB
SetPermanentCurrentInMicroAmp, 395	CFilterCoefficientsNet::s_FilterAttributesNet, 710
SetWaveLengthInNanometer, 396	Mcs::Usb, 65
CNF GenDeviceNet, 396	CompareFirmware
~CNF_GenDeviceNet, 396	CMcsUsbFactoryNet, 302
CNF_GenDeviceNet, 396	CompareTo
Set_Values, 397	HeadStageIDType, 695
COctoPotDeviceNet, 397	CompensateElectrodeOffset
BurnAdcOffset, 398	CWarnerUssingFunctionNet, 632
BurnDacOffset, 398	Connect
_aa	23.11000

OM	0 11 1 15 11 445
CMcsUsbNet, 322, 323	GetImplantResult, 415
CRFFunctionNet, 449 ConnectDevice	GetImplantState, 415
	GetOnOff, 415
CRadioControledDevicesNet, 444 ConnectedImp	GetPowerStrength, 416 GetRTC, 416
CPositionImpDeviceNet, 421	GetStateDebugData, 416
•	
ConnectImp	GetStateEventData, 417
CPositionImpDeviceNet, 421 ConnectSlave	RFFunction, 420 SetImplantCurrentSetpoint, 417
	SetPowerStrength, 417
CGilsonDeviceNet, 169 ControlState	SetRTC, 417
HeadStageIDTypeState, 698	SetN10, 417 SetStateDebugData, 419
CornerFrequency	SetStateEventData, 419
CFilterPropertyNet, 142	SwitchOnOff, 419
CornerFrequencymHz	CPositionImpDeviceNet, 420
CFilterPropertyNet, 142	!CPositionImpDeviceNet, 421
Count	~CPositionImpDeviceNet, 421
CMcsUsbListNet, 317	ConnectedImp, 421
CPatchServerDeviceNet, 404	ConnectImp, 421
CPatchServerDeviceNet, 405	CPositionImpDeviceNet, 421
Sensor, 405	GetDeviceList, 421
CPathIdentDeviceNet, 405	GetImpld, 422
~CPathIdentDeviceNet, 406	GetRFFrequency, 422
CPathIdentDeviceNet, 406	SetDeviceList, 422
Measure, 406	SetImpld, 422
Set_Values, 406	SetRFFrequency, 423
CPedoterDeviceNet, 406	CPPCDeviceNet, 423
!CPedoterDeviceNet, 407	CPPCDeviceNet, 423
~CPedoterDeviceNet, 407	McsBus, 424
CPedoterDeviceNet, 407	McsBus MotorControl, 424
GetCommand, 407	McsBus_Sensor, 424
SetCommand, 407	PPCFunction, 424
CPeristalticPumpDeviceNet, 408	CPPCFunctionNet, 424
~CPeristalticPumpDeviceNet, 408	!CPPCFunctionNet, 426
CPeristalticPumpDeviceNet, 408	~CPPCFunctionNet, 426
McsBus_MotorControl, 409	CPPCFunctionNet, 425
CPgaDeviceNet, 409	FirePressurePulse, 426
∼CPgaDeviceNet, 409	GetAnalogVoltage, 426
ApplyGains, 410	GetAnalogVoltageRange, 426
CPgaDeviceNet, 409	GetDigitalIn, 428
DefineAmplification, 410	GetPressureRange, 428
DefineFrequencyRange, 410	GetPumpModeType, 428
DefineNumAmplifications, 410	GetPumpSpeedUnit, 429
DefineNumFrequencyRanges, 410	GetSupplyVoltage, 429
GetAmplification, 410	GetValveActive, 429
GetFrequencyRange, 410	IsBusy, 429
GetGain, 411	LoadPressure, 431
GetNumAmplifications, 411	MeasureReservoir, 431
GetNumFrequencyRanges, 411	SetAnalogVoltageRange, 431
SetGain, 411	SetPressureOffset, 431
CPositionIIDeviceNet, 411	SetPressureRange, 431
!CPositionIIDeviceNet, 413	SetPumpModeType, 432
∼CPositionIIDeviceNet, 413	SetPumpSpeedUnit, 432
CPositionIIDeviceNet, 413	SetValveActive, 432
GetCoilCommunication, 413	CPPS_DeviceNet, 433
GetDebugData, 413	CPPS_DeviceNet, 433
GetEventData, 414	McsBus, 433
GetImplantCurrentSetpoint, 414	McsBus_MotorControl, 433

McsBus_Sensor, 433	SetFrequency, 445
PPS_Function, 433	StillConnected, 445
CPPS_FunctionNet, 434	CreateSideband
CPPS_FunctionNet, 434, 435	CStimulusFunctionNet, 576
GetAnalogVoltage, 435	CreateWirelessHeadstageSerialNumberString
GetAnalogVoltages, 435	CWirelessBaseFunctionNet, 680
GetBubbleState, 435	CRegionOfInterestRect
GetDigitalIn, 435	CCMOSMeaDeviceNet::CRegionOfInterestRect,
GetPumpCouple, 435	445
GetPumpEnableSpeedRatio, 435	CRetinaLedDeviceNet, 446
GetPumpFastOnOff, 435	~CRetinaLedDeviceNet, 447
GetPumpFastSpeed, 435	AddLoopEntry, 447
GetPumpFunctionSpeeds, 436	AddTableEntry, 447
GetPumpManualOnOff, 436	ClearTable, 447
GetPumpMaxSpeed, 436	CRetinaLedDeviceNet, 447
GetPumpModeType, 436	GetTablepointer, 447
GetPumpSpeedRatio, 436	SetLED, 447
GetPumpSpeedUnit, 436	SetLumi, 447
GetSupplyVoltage, 436	SetPersistency, 447
GetUseBubble, 436	SetRepeat, 448
SetAnalogVoltages, 436	SetTablepointer, 448
SetPumpCouple, 437	SetTrigger, 448
SetPumpEnableSpeedRatio, 437	CRFFunctionNet, 448
SetPumpFastOnOff, 437	!CRFFunctionNet, 449
SetPumpFastSpeed, 437	~CRFFunctionNet, 449
SetPumpFunctionSpeeds, 437	Connect, 449
SetPumpManualOnOff, 437	CRFFunctionNet, 449
SetPumpMaxSpeed, 437	GetAvailableDeviceList, 450
SetPumpModeType, 437	GetAvailableDeviceListEx, 450
SetPumpSpeedRatio, 438	GetAvailableStateList, 450
SetPumpSpeedUnit, 438	GetAvailableStateListEx, 450
SetUseBubble, 438	GetBaseFrequency, 451
CPPSDeviceNet, 438	GetConnectedDevice, 451
CPPSDeviceNet, 439	GetState, 451
CProgramPressureCurveNet, 439	GetTestMode, 451
!CProgramPressureCurveNet, 440	GetWorkingFrequency, 452
~CProgramPressureCurveNet, 439	SetBaseFrequency, 452
CProgramPressureCurveNet, 439	SetTestMode, 452
GetRepeats, 440	SetWorkingFrequency, 452
Program, 440	CRobo_FYIProgram_FunctionNet, 453
SetRepeats, 440	CRobo_FYIProgram_FunctionNet, 453
CPulseGeneratorFunctionNet, 440	GetLength, 453
!CPulseGeneratorFunctionNet, 441	GetState, 454
~CPulseGeneratorFunctionNet, 441	GetValve1, 454
CPulseGeneratorFunctionNet, 441	GetValve2, 454
GetModeSelect, 442	SetLength, 454
GetPeriod, 442	SetValve1, 454
GetPulseLength, 442	SetValve2, 454
SetModeSelect, 443	Start, 454
SetPeriod, 443	CRobo_FYITemp_FunctionNet, 454
SetPulseLength, 443	CRobo_FYITemp_FunctionNet, 455
CRadioControledDevicesNet, 443	GetICoeff, 455
ConnectDevice, 444	GetMaxPower, 455
CRadioControledDevicesNet, 444	GetPCoeff, 455
DisConnectDevice, 444	GetRegulatorOnOff, 455
GetDeviceNames, 444	GetSollTemp, 455
GetFrequency, 445	SetlCoeff, 456
HasRadioControl, 445	SetMaxPower, 456
·	*

CatDCaa# 4FC	CatDaciatamas C. 404
SetPoort, 456	GetResistanceC, 464
SetRegulatorOnOff, 456	GetResistanceV, 464
SetSollTemp, 456	GetScreen, 464
CRoboDacqNet, 456	Get Simulation, 464
CancelTableLoop, 459 CancelTableLoopAndStopTable, 459	GetUC, 464
·	GetUClamp, 465 GetUCOffset, 465
ClampAmpRestart, 459	
CRoboDacqNet, 459	GetUpdateDisplay, 465
DoRamp, 459 Emu GetCellCapacity, 460	GetUV, 465 GetUVOffset, 465
_ , , , ,	•
Emu_GetCellPotential, 460 Emu_GetCellResists, 460	GetXGain, 465
Emu_GetElectrodeResists, 460	RunTable, 465 SetAllDigout, 465
Emu_GetNoise, 460	SetCommand, 465
Emu_SetCellCapacity, 460	SetConfigurationBit, 465
Emu_SetCellPotential, 460	SetConfigurationBitAxc, 466
Emu_SetCellResists, 460	SetConfigurationBitAkc, 400 SetConfigurationBitBlu_Led, 466
Emu_SetElectrodeResists, 460	SetConfigurationBitBlu_LedToggleFast, 466
Emu_SetNoise, 460	SetConfigurationBitBlu_LedToggleSlow, 466
GetAllDigout, 461	SetConfigurationBitCC_Gen, 466
GetCapacityC, 461	SetConfigurationBitCV_Gen, 466
GetCapacityV, 461	SetConfigurationBitRC_Gen, 466
GetCapacityX, 461	SetConfigurationBitRed_Led, 466
GetClampAmpSerialNumber, 461	SetConfigurationBitRed_LedSaturation, 466
GetCommand, 461	SetConfigurationBitRed_LedToggleFast, 467
GetConfigurationBit, 461	SetConfigurationBitRed_LedToggleSlow, 467
GetConfigurationBitAxc, 461	SetConfigurationBitRelais, 467
GetConfigurationBitBlu_Led, 461	SetConfigurationBitRV_Gen, 467
GetConfigurationBitBlu_LedToggleFast, 461	SetConfigurationBitStream, 467
GetConfigurationBitBlu_LedToggleSlow, 461	SetConfigurationBitSupply, 467
GetConfigurationBitCC_Gen, 462	SetCrossTalkOffset, 467
GetConfigurationBitCV_Gen, 462	SetCrossTalkOptimum, 467
GetConfigurationBitRC_Gen, 462	SetDigout, 467
GetConfigurationBitRed Led, 462	SetDisplayText, 468
GetConfigurationBitRed_LedSaturation, 462	SetDownsampleFactor, 468
GetConfigurationBitRed LedToggleFast, 462	SetFilter, 468
GetConfigurationBitRed_LedToggleSlow, 462	SetFilterCoeffs, 468
GetConfigurationBitRelais, 462	SetIClamp, 468
GetConfigurationBitRV Gen, 462	SetICOffset, 468
GetConfigurationBits, 462	SetIGain, 468
GetConfigurationBitStream, 462	SetNoFilterCoeffs, 468
GetConfigurationBitSupply, 463	SetPGain, 468
GetCrossTalkOffset, 463	SetRecordingNumber, 469
GetCrossTalkOptimum, 463	SetScreen, 469
GetDigout, 463	SetSimulation, 469
GetDisplayText, 463	SetTriggerMaskValue, 469
GetDownsampleFactor, 463	SetUClamp, 469
GetFilter, 463	SetUCOffset, 469
GetFilterCoeffs, 463	SetUVOffset, 469
GetIC, 463	SetXGain, 469
GetlClamp, 463	StopTable, 469, 470
GetICOffset, 463	Table_Wait, 470
GetlGain, 464	TableDefBegin, 470
GetNIC_MS, 464	TableDefEnd, 470
GetNUC_MS, 464	TriggerMask_Default, 470
GetNUV_MS, 464	TriggerValue_MoveAbs, 470
GetPGain, 464	TriggerValue_StartQueue, 470
GetRecordingNumber, 464	UpdateDisplay, 470

VirtualDevice_ContinousDacq, 470	RoboError_CommandNotPossible, 484
VirtualDevice_TableRun, 470	RoboError_CommunicationTimeout, 484
CRoboDeviceNet, 471	RoboError_DacqNotReady, 484
∼CRoboDeviceNet, 474	RoboError_DLLMovementTimeout, 484
Axes_I, 482	RoboError_FindReferenceMethod, 484
Axes_X, 482	RoboError_GilsonCommandPending, 484
Axes_Y, 482	RoboError_GilsonTimeout, 484
Axes_Z, 483	RoboError_GilsonWrondID, 485
Axis_I, 483	RoboError_McsBus_UnknownCommand, 485
Axis X, 483	RoboError NoEndSwitch, 485
Axis_Y, 483	RoboError NoMoreData, 485
Axis_Z, 483	RoboError NoReference, 485
CancelPoolLoop, 474	RoboError_NoSpeedOrAcceleration, 485
CancelPoolLoopAndStopMovement, 474	RoboError_OverPressure, 485
CRoboDeviceNet, 474	RoboError_ParameterNotAllowed, 485
EnableQueue, 475	RoboError_PeristalticTimeout, 485
FindReference, 475	RoboError_Phase0OutOfRange, 485
GetAirpressure, 475	RoboError_PollLoopCanceled, 486
GetAirpressureLimit, 475	RoboError PollLoopCanceledAndStopMovement,
GetAirValve, 475	486
GetCurrentAirvalve, 475	
GetCurrentAirvalve, 475 GetCurrentAirvalveLimit, 475	RoboError_Pressure, 486 RoboError_RangeExceeded, 486
	_ -
GetCurrentPosition, 476	RoboError_StateChangeNotPossible, 486
GetErrorAirpressure, 476	RoboError_Timeout, 486
GetErrorCurrentAirvalve, 476	RoboError_UnknownCommand, 486
GetErrorVoltage12V, 476	RoboMainLowLevelCommand, 487
GetErrorVoltage5V, 476	RoboStatusEvent, 487
GetErrorVoltageAirvalve, 476	SetAirpressureLimit, 480
GetErrorVoltageRs485A, 476	SetAirValve, 480
GetErrorVoltageRs485B, 476	SetCurrentAirvalveLimit, 480
GetErrorVoltageValves, 477	SetCurrentAndAir, 480
GetInMovement, 477	SetInMovement, 480
GetMinPressure, 477	SetMinPressure, 481
GetMovementError, 477	SetVoltage12VLimit, 481
GetVoltage12V, 477	SetVoltage5VLimit, 481
GetVoltage12VLimit, 477	SetVoltageAirvalveLimit, 481
GetVoltage5V, 477	SetVoltageRs485ALimit, 481
GetVoltage5VLimit, 477	SetVoltageRs485BLimit, 481
GetVoltageAirvalve, 477	SetVoltageValvesLimit, 481
GetVoltageAirvalveLimit, 477	StartQueue, 481
GetVoltageRs485A, 478	StopMovement, 482
GetVoltageRs485ALimit, 478	WaitTimer, 482
GetVoltageRs485B, 478	CRoboDeviceNet::RoboMainLowLevelCommands, 702
GetVoltageRs485BLimit, 478	FindReferencePhase0, 703
GetVoltageValves, 478	GetAxisConfig, 703
GetVoltageValvesLimit, 478	GetHWConfig, 703
IsQueueEnabled, 478	GetHWRevision, 703
IsQueueStarted, 478	GetMaxNoPressure, 703
McsBus, 486	GetMaxNoPressureWaitTime, 703
McsBus_MotorControl, 486	GetMaxPressureWaitTime, 703
McsBus_XY, 483	GetMinNoPressureWaitTime, 703
McsBus_ZI, 483	GetMinPressure, 703
MoveAbs, 478, 479	GetMinPressureWaitTime, 704
NullCommand, 480	GetParameter, 704
RoboError_AnotherMaster, 483	GetPhases, 704
RoboError_Base, 484	GetSearchReferenceFastAccel, 704
RoboError_CannotEscapeEndSwitch, 484	GetSearchReferenceFastSpeed, 704
RoboError_CommandAlreadyInProgress, 484	GetSearchReferenceFineAccel, 704

GetSearchReferenceMethod, 705 GetSearchReferenceMoveOut, 705 GetSearchReferenceOffsetPos, 705 GetSearchReferenceOffsetPos, 705 GetUserParameter, 705 HasRef, 706 GetUserParameter, 705 HasRef, 706 SetAxisConfig, 706 SetAxisConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMessureWaitTime, 706 SetSetMinPressureWaitTime, 706 SetSetMinPressureWaitTime, 706 SetSetReferenceFineSpeed, 707 SetSearchReferenceFineAccel, 707 SetSearchReferenceFineAccel, 707 SetSearchReferenceFineAccel, 707 SetSearchReferenceFineSpeed, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetUnidDeviceNet, 487 ~CRoboFluidDeviceNet, 488 StopMovementZ, 498 StopMovementZ, 498 StopMovementZ, 498 StopMovementZ, 498	
GetSearchReferenceOffsetPos, 705 GetUserParameter, 705 HasRef, 706 SetAxisConfig, 706 SetAxisConfig, 706 SetHWConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativel, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetVelocityZ, 497 SetSpeedX, 497 SetVelocityZ, 497 SetVelocityZ, 497	
GetSearchReferenceOffsetPos, 705 GetUserParameter, 705 HasRef, 706 SetAxisConfig, 706 SetAxisConfig, 706 SetHWConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativel, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetSpeedX, 497 SetSpeedX, 497 SetVelocityZ, 497 SetVelocityZ, 497 SetSpeedX, 497 SetVelocityZ, 497 SetVelocityZ, 497	
GetUserParameter, 705 HasRef, 706 SetAxisConfig, 706 SetAxisConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedl, 497 SetSpeedNativel, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXy, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineSpeed, 707 SetSpeedXy, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityXy, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetDopNovementXy, 498 CRoboFluidDeviceNet, 487	
HasRef, 706 SetAxisConfig, 706 SetAxisConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSetMaxPressureWaitTime, 706 SetSetMinPressureWaitTime, 706 SetSetParameter, 707 SetSpeedNativeJ, 496 SetParameter, 707 SetSearchReferenceFastAccel, 707 SetSpeedNativeJ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMoveOut, 708 SetVelocityJ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 StopMovementI, 498 CRoboFluidDeviceNet, 487	
SetAxisConfig, 706 SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedl, 497 SetSpeedNativeI, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSpeedXy, 497 SetSearchReferenceFineSpeed, 707 SetSpeedXy, 497 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 708 SetVelocityZ, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 StopMovementI, 498 CRoboFluidDeviceNet, 487	
SetHWConfig, 706 SetHWRevision, 706 SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressure, 706 SetMinPressure, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedl, 497 SetParameter, 707 SetSpeedNativel, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXY, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMoveOut, 708 SetVelocityI, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 SetVelocityZ, 497 SetUserParameter, 708 StopMovementI, 498 CRoboFluidDeviceNet, 487	
SetHWRevision, 706 SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressure, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedl, 497 SetSpeedNativeWay, 496 SetSpeedNativeWay, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXy, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSpeedNativeZ, 497 SetSpeedXy, 497 SetSearchReferenceFineAccel, 707 SetSpeedXy, 497 SetSpeedZy, 497 SetSearchReferenceMethod, 707 SetSpeedZy, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityI, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 StopMovementXy, 498	
SetMaxNoPressure, 706 SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetSpeedI, 497 SetSpeedNativeI, 497 SetSpeedNativeI, 497 SetSpeedNativeXy, 498	
SetMaxNoPressureWaitTime, 706 SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressure, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetParameter, 707 SetSpeedNativeI, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXY, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMoveOut, 708 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 SetVelocityZ, 497 SetUserParameter, 708 StopMovementI, 498 CRoboFluidDeviceNet, 487	
SetMaxPressureWaitTime, 706 SetMinNoPressureWaitTime, 706 SetMinPressure, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetParameter, 707 SetSearchReferenceFastAccel, 707 SetSearchReferenceFastSpeed, 707 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityI, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetDevarentXY, 498	
SetMinNoPressureWaitTime, 706 SetMinPressure, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetSpeedNativeI, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXY, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedXY, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityI, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetOmorementXY, 498	
SetMinPressure, 706 SetMinPressureWaitTime, 706 SetSpeedI, 497 SetParameter, 707 SetSpeedNativeI, 497 SetSpeedNativeV, 497 SetSearchReferenceFastAccel, 707 SetSpeedNativeXY, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedXY, 497 SetSearchReferenceMethod, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityI, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetpMovementXY, 498	
SetMinPressureWaitTime, 706 SetSpeedI, 497 SetParameter, 707 SetSearchReferenceFastAccel, 707 SetSearchReferenceFastSpeed, 707 SetSpeedNativeXY, 497 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetVelocityXY, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetMovementXY, 498	
SetParameter, 707 SetSearchReferenceFastAccel, 707 SetSearchReferenceFastSpeed, 707 SetSearchReferenceFastSpeed, 707 SetSpeedNativeZ, 497 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetSpeedNativeZ, 497 SetSpeedNativeZ, 497 SetSpeedNativeZ, 497 SetSpeedNativeZ, 497 SetSpeedNativeZ, 497 SetVelocityI, 497 SetVelocityJ, 497 SetVelocityZ, 497 SetOpMovementI, 498 CRoboFluidDeviceNet, 487	
SetSearchReferenceFastAccel, 707 SetSearchReferenceFastSpeed, 707 SetSearchReferenceFineAccel, 707 SetSearchReferenceFineAccel, 707 SetSpeedXY, 497 SetSearchReferenceFineSpeed, 707 SetSpeedZ, 497 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetSpeedXY, 497 SetVelocityI, 497 SetVelocityZ, 497 SetVelocityZ, 497 SetOpMovementI, 498 StopMovementXY, 498	
SetSearchReferenceFastSpeed, 707 SetSearchReferenceFineAccel, 707 SetSearchReferenceFineSpeed, 707 SetSearchReferenceFineSpeed, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityX, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetSpeedXY, 497 SetSpeedXY, 497 SetVelocityI, 497 SetVelocityZ, 497 SetVelocityZ, 497 StopMovementI, 498 StopMovementXY, 498	
SetSearchReferenceFineAccel, 707 SetSearchReferenceFineSpeed, 707 SetSearchReferenceFineSpeed, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceMoveOut, 708 SetVelocityXY, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 StopMovementI, 498 CRoboFluidDeviceNet, 487 SetVelocityX, 498	
SetSearchReferenceFineSpeed, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetVelocityXY, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 StopMovementI, 498 CRoboFluidDeviceNet, 487 SetVelocityZ, 497 StopMovementXY, 498	
SetSearchReferenceMethod, 707 SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetVelocityXY, 497 SetSearchReferenceOffsetPos, 708 SetVelocityZ, 497 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetVelocityZ, 497 StopMovementI, 498 StopMovementXY, 498	
SetSearchReferenceMoveOut, 708 SetSearchReferenceOffsetPos, 708 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetVelocityXY, 497 SetVelocityZ, 497 StopMovementI, 498 StopMovementXY, 498	
SetSearchReferenceOffsetPos, 708 SetUserParameter, 708 CRoboFluidDeviceNet, 487 SetVelocityZ, 497 StopMovementI, 498 StopMovementXY, 498	
SetUserParameter, 708 StopMovementI, 498 CRoboFluidDeviceNet, 487 StopMovementXY, 498	
CRoboFluidDeviceNet, 487 StopMovementXY, 498	
·	
CloseAllValves, 488 CRoboStatorDeviceNet::RoboMainStatorLowLevelComma	ands.
CRoboFluidDeviceNet, 488 709	,
GetPumpSpeed, 488 FindReferencePhase0XY, 709	
GetSingleValve, 488 CSafeISDeviceNet, 498	
GetValve, 488 ~CSafeISDeviceNet, 499	
IsPumpMotorOn, 489 CSafeISDeviceNet, 499	
m_pMcsBus_MotorControlNet, 490 DacqDevice, 501	
m_pRoboFluidDevice, 490 FluidControlDevice, 501	
McsBus_MotorControl, 490 RoboDevice, 501	
PumpOff, 489 SetAdcChannels, 500	
PumpOn, 489 SetAdcSamplePos, 500	
SetPumpSpeed, 489 SetDacMode, 500	
SetSingleValve, 489 SetDacPeriode, 500	
SetValve, 489 SetDacPulseform, 500	
CRobolnjectDeviceNet, 490 SetSwitches, 501	
CRobolnjectDeviceNet, 491 csCapacityTestDischarge	
CRoboocyte2DeviceNet, 491 Mcs::Usb, 69	
CRoboocyte2DeviceNet, 492 csCapacityTestPrecharge	
GetAxisLED, 492 Mcs::Usb, 69	
GetGilsonDevice, 492 csCharge	
GetMcsBus_Extension, 492 Mcs::Usb, 69	
GetRoboDacq, 492 CSCUDacqGroupChannelSelectionNet, 501	
GetRoboFluidDevice, 492 CSCUDacqGroupChannelSelectionNet, 502	
SetAxisLED, 492 CSCUFunctionNet, 502	
CRoboStatorDeviceNet, 492 !CSCUFunctionNet, 505	
CRoboStatorDeviceNet, 494 ~CSCUFunctionNet, 505	
FindReferencel, 494 AutomaticAnalogOut, 505	
FindReferenceXY, 494 CSCUFunctionNet, 505	
FindReferenceZ, 494 EnableAnalogOut, 505	
GetCurrentPositionI, 494 GetAnalogOutADCRange, 506	
GetCurrentPositionXY, 494 GetAnalogOutChannels, 506	
GetCurrentPositionZ, 495 GetAnalogOutDACRange, 506	

GetAvailableHeadstages, 506	ReceiveString, 520
GetAvailableHeadstagesEvent, 519	Send, 520
GetFilterProperties, 506	csError
GetFilterProperty, 507	Mcs::Usb, 69
GetHeadstageAdcBits, 507	csldleChargeFinished
GetHeadstageAdcRangeInMicroVolt, 507	Mcs::Usb, 69
GetHeadstageDacBits, 508	csldleNoBattery
GetHeadstageDacCurrentRangeInMicroAmpere,	Mcs::Usb, 69
508	csRefreshBattery
GetHeadstageDacCurrentResolutionInNanoAm-	Mcs::Usb, 69
pere, 508	CStg200xBasicNet, 520
GetHeadstageDacVoltageRangeInMilliVolt, 509	~CStg200xBasicNet, 525
GetHeadstageDacVoltageResolutionInMicroVolt,	GetAnalogRanges, 526
509	GetAnalogResolution, 526
GetHeadstageFrameCyclesToCompare, 509	GetAutocalibrationDisabled, 526
GetHeadstageGainInPermille, 510	GetAvailableMemory, 527
GetHeadstageID, 510	GetBlankingEnable, 527
GetHeadstageLinkSpeed, 510	GetCanDisableCalibration, 528
GetHeadstageNumberOfAnalogChannels, 511	GetCanDisableCurrentswitch, 528
GetHeadstageNumberOfStimulationChannels, 511	GetCurrentRangeInNanoAmp, 528
GetHeadstagePowerStateAtStart, 511	GetCurrentRangeInNanoAmpByIndex, 528
GetHeadstageSamplerate, 512	GetCurrentRangeListInNanoAmp, 529
GetHeadstageSerialNumber, 512	GetCurrentRangeSelectedIndex, 529
GetMaxNumberOfHeadstages, 512	GetCurrentResolutionInNanoAmp, 529
GetMaxStimulusChannelsPerHeadstage, 512	GetCurrentResolutionInNanoAmpByIndex, 530
GetReferenceElectrodeMode, 513	GetDacAmplificationFactor, 530
GetReferenceElectrodeSwitchState, 513	GetDACResolution, 530
HasAnalogOut, 513	GetDiginValue, 530
HasGalvanicIsolation, 513	GetDigntValue, 530 GetDigoutMode, 531
HasHSPowerSwitch, 514	GetDigoutWide, 531
IsAnalogOutEnabled, 514	GetElectrodeDacMux, 531
IsAutomaticAnalogOut, 514	GetElectrodeEnable, 532
IsHeadstageAvailable, 514	GetElectrodeEnable, 532
IsHeadstageAvailableEvent, 519	GetEnableAmplifierProtectionSwitch, 534
IsHSPowered, 515	GetExternalElectrodeEnable, 535
IsInDacqLegacyMode, 515	GetEAAmplification, 535
OnGetAvailableHeadstages, 515	GetHasLowCurrentswitchLatency, 535
OnlsHeadstageAvailable, 515	GetHeadstage, 535
PowerHS, 515	GetListmodeIndexRange, 536
SetAnalogOutADCRange, 516 SetAnalogOutChannels, 516	GetListmodeTriggerSource, 536
•	GetNumberOfAnalogChannels, 536
SetAnalogOutDACRange, 516	GetNumberOfCurrentRangeIndexes, 537
Set DacqLegacyMode, 516	GetNumberOfHWDACPaths, 537
SetHeadstageFrameCyclesToComparePermanent,	GetNumberOfStimulationElectrodes, 537
517	GetNumberOfStimulationSourcesPerElectrode,
SetHeadstageLinkSpeedPermanent, 517	537
SetHeadstageNumberOfAnalogChannelsPerma-	GetNumberOfSyncoutChannels, 537
nent, 517	GetNumberOfTriggerInputs, 537
SetHeadstagePowerStateAtStart, 517	GetNumberOfVoltageRangeIndexes, 538
SetHeadstageSampleratePermanent, 518	GetOutputRate, 538
SetReferenceElectrodeMode, 518	GetStgProgramInfo, 538
SetReferenceElectrodeSwitchState, 518	GetStgVersionInfo, 539
csDischarge	GetSyncoutMap, 539
Mcs::Usb, 69	GetTimeResolutionInNanoSeconds, 539
CSerialPortNet, 519	GetTotalMemory, 540
CSerialPortNet, 519	GetTriggerSource, 540
GetBytesAvailable, 520	GetVoltageRangeInMicroVolt, 540
Receive, 520	GetVoltageRangeInMicroVoltBvIndex, 540

GetVoltageRangeListInMicroVolt, 541	GetModuleTemp, 569
GetVoltageRangeSelectedIndex, 541	MwPollStatusEvent, 573
GetVoltageResolutionInMicroVolt, 541	PrepareAndAppendData, 569
GetVoltageResolutionInMicroVoltByIndex, 541	PrepareAndSendData, 570
HasDigitalPort, 542	QueryTriggerstatus, 571
HasFilterAmp, 542	SegmentDefine, 571
LegacyHasOutputMap, 542	SegmentSelect, 571
ListModeSendStart, 542	SegmentStart, 572
ListModeSendStop, 543	SetOutputMap, 572
SendStart, 543	Stg200xPollStatusEvent, 573
SendStop, 543	CStimulusFunctionNet, 573
SetAutocalibrationDisabled, 544	ClearChannel_PrepareAndSendData, 575
	_ ·
SetBlankingEnable, 544, 545	ClearChannelData, 575
SetCurrentBonge SelectedIndov, 546	ClearMultiplexedData, 575
SetCurrentRangeSelectedIndex, 546	ClearSyncData, 575
SetDacAmplificationFactor, 546	CreateSideband, 576
SetDigoutMode, 546	CStimulusFunctionNet, 575
SetDigoutValue, 546	ForceStatusEvent, 577
SetElectrodeDacMux, 547, 548	GetAvailableMemory, 577
SetElectrodeEnable, 549, 550	GetCurrentRangeInNanoAmp, 577
SetElectrodeMode, 551, 552	GetCurrentResolutionInNanoAmp, 577
SetEnableAmplifierProtectionSwitch, 553, 554	GetDACResolution, 578
SetExternalElectrodeEnable, 554, 555	GetMultiplexedDataChannelsInBlock, 578
SetFAAmplification, 555	GetNumberOfAnalogChannels, 578
SetHeadstage, 555	GetTotalMemory, 578
SetListmodeIndexRange, 555	GetVoltageRangeInMicroVolt, 578
SetListmodeTriggerSource, 556	GetVoltageResolutionInMicroVolt, 579
SetMeasurementMode, 557	PollStatusEvent, 584
SetOutputRate, 557	PrepareAndAppendData, 579
SetStgProgramInfo, 557	PrepareAndSendData, 580
SetSyncoutMap, 557	PrepareData, 581
SetTriggerSource, 558	SendMultiplexedData, 581
SetVoltageMode, 558	SendPreparedData, 581
SetVoltageRangeSelectedIndex, 558	SendStart, 581
CStg200xDownloadBasicNet, 559	SendStop, 582
ClearChannelData, 561	SetupTrigger, 582
ClearSyncData, 561	SetupTriggerSingle, 583
DisableAutoReset, 561	StartPoll, 583
EnableAutoReset, 561	StopPoll, 583
ForceStatusEvent, 561	CStimulusFunctionNet::SidebandData, 711
GetMemoryUsageDAC, 561	!SidebandData, 711
GetMemoryUsageSyncout, 562	~SidebandData, 711
GetSweepCount, 562	Duration, 711
GetTrigger, 562	Sideband, 711
ResetStatus, 563	SidebandData, 711
SendChannelData, 563	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData
SendSyncData, 564	712
SetupRetriggerMode, 564	!StimulusDeviceDataAndUnrolledData, 713
	~StimulusDeviceDataAndUnrolledData, 713
SetupTrigger, 565	
SetupTriggerSingle, 565	DeviceData, 713
Stimulus, 566	DeviceDataLength, 713
CStg200xDownloadNet, 566	StimulusDeviceDataAndUnrolledData, 713
~CStg200xDownloadNet, 568	UnrolledAmplitude, 713
ClearChannel_PrepareAndSendData, 568	UnrolledDuration, 713
CStg200xDownloadNet, 567, 568	UnrolledSync, 713
DisableMultiFileMode, 569	CSw2to64DeviceNet, 584
EnableMultiFileMode, 569	~CSw2to64DeviceNet, 585
GetModuleCurrent, 569	CSw2to64DeviceNet, 585

GetChannel, 585	GetSetpointMax, 595
GetChannels, 585	GetSetpointMin, 595
GetNumber, 585	GetThermocoupleCalibration, 595
SetChannel, 585	GetThermocoupleNanovoltPerKelvin, 595
SetChannels, 586	GetThermocoupleReferenceTemp, 596
CTcxDeviceNet, 586	GetThermocoupleTemp, 596
∼CTcxDeviceNet, 589	GetThermocoupleTempAbs, 596
CalibrateThermocouple, 589	GetUnit, 596
CTcxDeviceNet, 588	GetUOut, 596
FactoryReset, 589	GetValue, 596
GetBoardTemp, 589	GetValueHires, 597
GetCalibration, 589	GetVolti, 597
GetCalibrationDecp, 589	SetCalibration, 597
GetCalibrationMax, 589	SetD, 597
GetCalibrationMin, 589	SetDevice, 597
GetCurrent, 589	SetDeviceType, 597
GetD, 590	SetDevname, 597
,	SetEnableHeaterLimit, 597
GetDecp, 590	•
GetDevice, 590	SetEnableThermocouple, 598
GetDeviceType, 590	SetHeaterLimit, 598
GetDevname, 590	Setl, 598
GetDMax, 590	SetMaxHeaterPowerMultiwell, 598
GetDMin, 590	SetMaxP, 598
GetDuty, 590	SetOnOff, 598
GetEnableHeaterLimit, 591	SetP, 599
GetEnableThermocouple, 591	SetSensorType, 599
GetHasThermocouple, 591	SetSetpoint, 599
GetHeaterLimit, 591	SetThermocoupleNanovoltPerKelvin, 599
GetHeaterTemp, 591	CTEERFunctionNet, 599
Getl, 591	!CTEERFunctionNet, 602
GetIDecp, 591	\sim CTEERFunctionNet, 602
GetlMax, 592	CancelInternalCalibration, 602
GetlMin, 592	CTEERFunctionNet, 601
GetlOut, 592	GetAdapterCode, 602
GetMaxHeaterPowerMultiwell, 592	GetAdcOffsetU1, 602
GetMaxP, 592	GetAdcOffsetU2, 602
GetMaxpDecp, 592	GetAmplitude_nA, 602
GetMaxpMax, 592	GetBytesPerSample, 603
GetMaxpMin, 593	GetClampMode, 603
GetNumControlChannels, 593	GetControllerParams, 603
GetNumDevices, 593	GetCurrentEnable, 603
GetNumMeasureChannels, 593	GetDacZero, 604
GetOnOff, 593	GetFrameErrorCounter, 604
GetP, 593	GetLiquidResistance, 604
GetPDecp, 593	GetMaxChunkSize Byte, 604
GetPMax, 593	GetNumberOfAvailableSamples, 604
	•
GetPOirt 503	GetPeriod_us, 605
GetPOut, 593	GetRotaryPositionCode, 605
GetPwrOut, 594	GetSampleBufferChunk, 605
GetPwrSet, 594	GetSampleRate, 605
GetRes1, 594	GetSampleVoltageBuffer_uV, 606
GetRes2, 594	GetScaleFactorU1, 606
GetResS, 594	GetScaleFactorU2, 606
GetResX, 594	GetUptimeSeconds, 606
GetROut, 594	GetWaveform, 606
GetSensorType, 595	IsInternalCalibrationFinished, 607
GetSetpoint, 595	IsSamplingFinished, 607
GetSetpointDecp, 595	SetAmplitude_nA, 607

SetBufferIndex, 607	GetHeadstageSamplingActive, 618
SetClampMode, 608	GetMultiHeadstageMode, 618
SetControllerParams, 608	GetPicFirmwareType, 618
SetCurrentEnable, 608	GetSelectedChannels, 619
SetExternalLED, 608	GetSelectedHeadstageState, 619
SetLiquidResistance, 609	GetStimulusParametersCache, 619
SetPeriod_us, 609	GetStimulusParametersFromSelectedHS, 619
SetWaveform, 609	GetStiumlusParameters, 619
StartInternalCalibration, 609	GetUserDefinedName, 619
StartSampling, 609	GetUserDefinedNameCache, 619
StopSampling, 610	GetUserDefinedNameFromSelectedHS, 620
CTEERMachineDeviceNet, 610	PulseGenerator, 621
~CTEERMachineDeviceNet, 610	SelectHeadstage, 620
CTEERMachineDeviceNet, 610	SetAccelGyroDesiredRate, 620
TEERFunctionNet, 611	SetAccelGyroEnabled, 620
CurrentClamp	SetAccelRange, 620
Mcs::Usb, 90	SetAnalogOutChannel, 620
CurrentMeasure	SetAnalogOutFilter, 620
Mcs::Usb, 57	SetAudioChannels, 620
CurrentRangeInNanoAmp	SetDacRange, 621
- ·	-
W2100_StimulusParametersNet, 715	SetUpedetage On Off 631
CurrentResolutionInNanoAmp	SetHeadstageOnOff, 621
W2100_StimulusParametersNet, 715	SetHeadstageSamplingActive, 621
CUsbDeviceConfigurationFunctionNet, 611	SetHeadstageToSleep, 621
!CUsbDeviceConfigurationFunctionNet, 612	SetMultiHeadstageMode, 621
~CUsbDeviceConfigurationFunctionNet, 612	SetSelectedChannels, 621
CUsbDeviceConfigurationFunctionNet, 611	Stimulator, 622
	CW2100_FunctionNet::AudioChannelsNet, 96
SetDeviceName, 612	amplification, 96
CUsbExceptionNet, 612	channel, 96
CUsbExceptionNet, 613	dacqgroup, 96
	CW2100_StimulatorFunctionNet, 622
CutoffOrCenterFrequency	BOOST_BIT, 628
CCreateFilterNet, 118	ClearChannelData, 623
CVoltageRangeInfoNet	CW2100_StimulatorFunctionNet, 623
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	, GetBoostAlwaysOnMode, 624
614	GetBoostPreTime, 624
CW2100_FunctionNet, 614	GetCurrentRangeInNanoAmp, 624
ClearStimulusParametersCache, 616	GetCurrentResolutionInNanoAmp, 624
ClearUserDefinedNameCache, 616	GetDACResolution, 624
CW2100_FunctionNet, 616	GetDigitalStimulatorTrigger, 624
DeselectAllHeadstages, 617	GetDigitalStimulatorTriggerSlope, 625
DeselectHeadstage, 617	GetNumberOfAnalogChannels, 625
GetAccelGyroCurrentRate, 617	GetNumberOfSyncoutChannels, 625
GetAccelGyroDesiredRate, 617	GetNumberOfTriggerInputs, 625
GetAccelGyroEnabled, 617	GetStimulationPatternMemory, 625
GetAccelRange, 617	GetTimeResolutionInNanoSeconds, 625
GetAnalogOutChannel, 617	GetTimeSlot, 625
GetAnalogOutFilter, 617	GetVoltageRangeInMicroVolt, 625
GetAudioChannels, 617	GetVoltageResolutionInMicroVolt, 626
GetAvailableHeadstages, 617	GND_SWITCH_BIT, 628
GetBatteryState, 618	PollStatusEvent, 628
GetDacRange, 618	PrepareData, 626
GetBackange, 618 GetFilterProperties, 618	PrepareDataSync, 626
•	SelectTimeSlot, 626
GetFilterProperty, 618	
GetFPGAFirmwareType, 618	SendPreparedData, 626
GetGyroRange, 618	SendStart, 627
GetHeadstageOnOff, 618	SendStop, 627

SetDigitalStimulatorTrigger, 627	WaitForChamber, 646
SetDigitalStimulatorTriggerSlope, 627	CWarnerValveControllerDeviceNet, 646
StartPoll, 627	!CWarnerValveControllerDeviceNet, 651
StopPoll, 627	~CWarnerValveControllerDeviceNet, 651
SYNC_BIT0, 628	ClearTableName, 651
SYNC BIT1, 628	ClearValveTable, 651
CW2100DacqGroupChannelSelectionNet, 628	CWarnerValveControllerDeviceNet, 651
CW2100DacqGroupChannelSelectionNet, 628	GetActiveRunningTableNumber, 651
CWarnerUssingDeviceNet, 629	GetActiveRunningTableNumberEvent, 668
!CWarnerUssingDeviceNet, 630	GetAnalogThresholdHigh, 651
~CWarnerUssingDeviceNet, 629	GetAnalogThresholdHighEvent, 668
CWarnerUssingDeviceNet, 629	GetAnalogThresholdLow, 652
WarnerUssingFunction, 630	GetAnalogThresholdLowEvent, 668
CWarnerUssingFunctionNet, 630	GetAnalogVoltage, 652
!CWarnerUssingFunctionNet, 632	GetAnalogVoltageEvent, 668
~CWarnerUssingFunctionNet, 632	GetCurrentEditTableNumber, 652
CompensateElectrodeOffset, 632	GetCurrentNumberOfValves, 652
CWarnerUssingFunctionNet, 632	GetCurrentNumberOfValvesEvent, 669
GetAvailableChambers, 633	GetDigitalOutPortValve, 653
GetChannelsCountOfChamber, 633	GetDigitalOutPortValveEvent, 669
GetClampMode, 633	GetDigitalPortDirection, 653
GetComplianceVoltageThreshold, 633	GetDigitalPortDirectionEvent, 669
GetDacPampsPerDigitHighCurrentRange, 634	GetDisplayMode, 653
GetDacPampsPerDigitLowCurrentRange, 634	GetDisplayModeEvent, 669
GetDacZero, 634	GetTableName, 653
GetHighCurrentRange, 635	GetTableNamebyIndex, 654
GetIdleModeOffset, 635	GetTableNamebyIndexEvent, 669
GetLiquidResistance, 635	GetTotalNumberOfDigitalPorts, 654
GetLowCurrentRange, 636	GetTotalNumberOfTables, 654
GetNumberOfAvailableChambers, 636	GetTotalNumberOfValves, 654
GetNumberOfHardwareSlotsForChambers, 636	GetTotalTableSize, 655
GetU1Offset, 636	GetValveActive, 655
GetU1Reference, 638	GetValveActiveEvent, 669
GetU2Offset, 638	GetValveBoardRevision, 655
GetU2Reference, 638	GetValveBoardRevisionEvent, 669
GetUnitDescription, 639	GetValveBoardRevisionString, 655
GetUnitExponent, 639	GetValveCurrent, 655
GetUnitName, 639	GetValveDigitalInPort, 656
GetUnitsPerDigit, 640	GetValveDigitalInPortEvent, 670
GetUptimeSeconds, 640	GetValveLedOn, 656
GetVoltageClampControllerParam_D, 640	GetValveLedOnEvent, 670
GetVoltageClampControllerParam_I, 641	GetValveManualGroup, 656
GetVoltageClampControllerParam_P, 641	GetValveManualGroupEvent, 670
IsChamberAvailable, 641	GetValveManualState, 657
IsHighCurrentMode, 642	GetValveManualStateEvent, 670
IsInternalCalibrationFinished, 642	GetValveMode, 657
IsPulseEnabled, 642	GetValveModeEvent, 670
SetClampMode, 643	GetValvesActiveMap, 657
SetEnablePulse, 643	GetValvesManualStateMap, 657
SetHighCurrentMode, 643	GetValveTableEntry, 658
SetIdleModeOffset, 644	IsDigitalOutPortInverted, 658
SetLiquidResistance, 644	IsDigitalOutPortInvertedEvent, 670
SetLowCurrentMode, 644	IsValveDigitalInInverted, 658
SetPulse, 644	IsValveDigitaIInInvertedEvent, 670
SetVoltageClampControllerParam_D, 645	IsValveOpen, 659
SetVoltageClampControllerParam_I, 645	IsValveOpenEvent, 671
SetVoltageClampControllerParam_P, 645	IsValveOpenInAnalogMode, 659
WaitForAllChambers, 646	IsValveOpenInAnalogModeEvent, 671

	1.V.I. O I. D.: '! IM . I 050	0.10 0.70
	IsValveOpenInDigitalMode, 659	GetSync, 673
	IsValveOpenInDigitalModeEvent, 671	SetADC, 673
	LoadValveTable, 659	SetIO, 673
	OnGetActiveRunningTableNumber, 660	SetIODirection, 674
	OnGetAnalogThresholdHigh, 660	SetTrigger, 674
	OnGetAnalogThresholdLow, 660	SetTriggerSyncDirection, 674
	OnGetAnalogVoltage, 660	CWClassicFunctionNet, 674
	OnGetCurrentNumberOfValves, 660	CWClassicFunctionNet, 675, 676
	OnGetDigitalOutPortValve, 660	GetFilterParametersHeadstage, 676
	OnGetDigitalPortDirection, 660	GetHasChecksum, 676
	OnGetDisplayMode, 660	GetHasRedLedHeadstage, 676
	OnGetTableNamebyIndex, 660	GetHeadstageOnOff, 676
	OnGetValveActive, 661	GetResetFilter, 676
	OnGetValveBoardRevision, 661	GetRFConnectionStatus, 676
	OnGetValveDigitalInPort, 661	GetRFFrequencyHeadstage, 676
	OnGetValveLedOn, 661	GetRFFrequencyReceiver, 676
	OnGetValveManualGroup, 661	GetRFPower, 677
	OnGetValveManualState, 661	GetScanHeadstagesResult, 677
	OnGetValveMode, 661	GetSelectedHeadstage, 677
	OnIsDigitalOutPortInverted, 661	GetSerialNumberHeadstage, 677
	OnIsValveDigitalInInverted, 662	GetWPADebugMode, 677
	OnIsValveOpen, 662	GetWPAType, 677
	OnIsValveOpenInAnalogMode, 662	ResetChannelmap, 677
	OnIsValveOpenInDigitalMode, 662	ScanForHeadstages, 677
	OnTableEntryChanged, 662	SetChannelmap, 677
	SetActiveRunningTableNumber, 662	SetFilterParametersHeadstage, 677
	SetAnalogThresholdHigh, 662	SetHasChecksum, 678
	SetAnalogThresholdLow, 663	SetHeadstageOnOff, 678
	SetCurrentEditTableNumber, 663	SetHWSelectedChannels, 678
	SetDefault, 663	SetResetFilter, 678
	SetDigitalOutPortInvert, 663	SetRFFrequencyHeadstage, 678
	SetDigitalOutPortValve, 664	SetRFFrequencyReceiver, 678
	SetDigitalPortDirection, 664	SetRFFrequencyReceiverEeprom, 678
	SetDisplayMode, 664	SetRFLostBehaviour, 678
	SetTableName, 664	SetRFPower, 679
	SetTableStep, 665	SetSelectedHeadstage, 679
	SetTableStepAll, 665	SetSerialNumberHeadstage, 679
	SetValveActive, 665	SetWPADebugMode, 679
	SetValveCurrent, 665	SetWPAType, 679
	SetValveDigitalInInvert, 666	CWirelessBaseFunctionNet, 679
	SetValveDigitalInPort, 666	CreateWirelessHeadstageSerialNumberString,
	SetValveLedOn, 666	680
	SetValveManualGroup, 666	CWirelessBaseFunctionNet, 680
	SetValveManualState, 667	CyclePort
	SetValveMode, 667	CMcsUsbNet, 324
	SetValvesActiveMap, 667	Cypress
	SetValvesManualStateMap, 667	Mcs::Usb, 91
	SetValveTableEntry, 667	Cypress FX1
	StoreValveTable, 668	Mcs::Usb, 77
	TableEntryChangedEvent, 671	Cypress_FX2
CW	arnerValveControllerDeviceTesterFunctionNet, 671	Mcs::Usb, 77
Ovvo	!CWarnerValveControllerDeviceTesterFunctionNet,	Cypress_FX3
	673	Mcs::Usb, 77
	~CWarnerValveControllerDeviceTesterFunctionNet,	1VIU3U3D, 11
	672	DAC1Channel
	CWarnerValveControllerDeviceTesterFunctionNet,	Mcs::Usb, 60
	672	DAC2Channel
	GetIO, 673	Mcs::Usb, 60
	Getto, 0/3	DAC3Channel

Mcs::Usb, 60	DEST_TARGET5
DAC4Channel	Mcs::Usb, 56
Mcs::Usb, 60	DEST_TARGET6
DACQ1 DigitalGroup	Mcs::Usb, 56
Mcs::Usb, 61	DEST_TARGET7
DacqDevice	Mcs::Usb, 56
CSafeISDeviceNet, 501	DEST_TARGET8
dacqgroup	Mcs::Usb, 56
CW2100_FunctionNet::AudioChannelsNet, 96	DEST_TARGET9
DacqGroupChannelEnumNet	Mcs::Usb, 56
Mcs::Usb, 60	DEST_TARGET_MASK
DacqMeaGroupTypeEnumNet	Mcs::Usb, 56
Mcs::Usb, 61	DetectChipType
DacqTrigger	CCMOSMea_FunctionNet, 104
Mcs::Usb, 64	DEVICE_NOT_FOUND
DACResolution	Mcs::Usb, 65
W2100_StimulusParametersNet, 715	DeviceArrival
DataModeEnumNet	CMcsUsbListNet, 317
Mcs::Usb, 61	Device Data CStimulus Function Netw Stimulus Davis a Data And Haralled Data
DataState LlandState DT:::::aState COO	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
HeadStageIDTypeState, 698 DeepCopy	713
	DeviceDataLength CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
CCMOSMeaDeviceNet::CRegionOfInterestRect, 445	713
DefineAmplification	DeviceEnumNet
CPgaDeviceNet, 410	Mcs::Usb, 61
DefineFrequencyRange	DeviceHasNoHeadstage
CPgaDeviceNet, 410	Mcs::Usb, 67, 68, 82
DefineNumAmplifications	DeviceId
CPgaDeviceNet, 410	CMcsUsbListEntryNet, 314
DefineNumFrequencyRanges	DeviceIdNet, 680
CPgaDeviceNet, 410	BcdDevice, 681
DeselectAllHeadstages	BusType, 681
CW2100 FunctionNet, 617	DeviceIdNet, 681
DeselectHeadstage	IdProduct, 681
CW2100_FunctionNet, 617	IdVendor, 681
DEST_FX3_TARGET_MASK	operator=, 681
Mcs::Usb, 56	DeviceName
DEST_TARGET1	CMcsUsbListEntryNet, 314
Mcs::Usb, 56	DeviceNotConnected
DEST_TARGET10	Mcs::Usb, 66, 67, 81
	DeviceRemoval
DEST_TARGET11	CMcsUsbListNet, 317
	DeviceRunStatus
DEST_TARGET12	Mcs::Usb, 57, 71, 74, 75, 83, 86, 93
	DigDataFromReceiver
DEST_TARGET13	Mcs::Usb, 93
Mcs::Usb, 56	Digital
DEST_TARGET14	Mcs::Usb, 76, 94
Mcs::Usb, 56	DigitalData
DEST_TARGET15	Mcs::Usb, 57, 71, 74, 75, 83, 86, 93
Mcs::Usb, 56	DigitalDatastreamEnableEnumNet
DEST_TARGET2	Mcs::Usb, 63
Mcs::Usb, 56	DigitalGroup
DEST_TARGET3	Mcs::Usb, 61
Mcs::Usb, 56	DigitalIn
DEST_TARGET4	Mcs::Usb, 57, 63, 71, 73, 74, 83, 86, 93
Mcs::Usb, 56	DigitalInOfOutPort

Mcs::Usb, 57, 71, 73, 74, 83, 86, 93	Disconnect
DigitalInPort	CMcsUsbNet, 324
Mcs::Usb, 59	DisConnectDevice
DigitalInReserverd	CRadioControledDevicesNet, 444
Mcs::Usb, 63	DongleS
DigitalMux	Mcs::Usb, 78
Mcs::Usb, 59	DoRamp
DigitalOut	CRoboDacqNet, 459
Mcs::Usb, 63	Dotriapot
DigitalOutReg	Mcs::Usb, 77
Mcs::Usb, 59	DoubleToInt
DigitalOutReserved	Mcs::Usb, 66
Mcs::Usb, 63	DownloadFirmware
DigitalOutStimulator	CMcsUsbFactoryNet, 302
Mcs::Usb, 57, 71, 74, 75, 83, 86, 93	DownloadOnly
DigitalPulse	Mcs::Usb, 85
Mcs::Usb, 57, 71, 73, 74, 83, 86, 93	DriverVersionNet, 684
DigitalReg	\sim DriverVersionNet, 685
Mcs::Usb, 59	DriverVersionNet, 685
DigitalSource	DriverVersionNet::FormatVersion, 685
DigitalSource< digitalsourceenum >, 682	GetDestinationCode, 686
DigitalSource< digitalsourceenum >, 682	GetDestinationName, 687
DigitalSource, 682	GetMajor, 687
MaxBitNumber, 682	GetMinor, 688
MaxBitNumberStatic, 682	GetNumEntries, 688
size, 682	GetSerialNumber, 688
Source, 683	GetStatus, 689
DigitalSourceGeneral, 683	GetVersionInt, 689
DigitalSourceGeneral, 683	GetVersionString, 690
MaxBitNumber, 683, 684	DriverVersionNet::FormatVersion
size, 684	DriverVersionNet, 685
Source, 684	DSP
DigitalStimulatorTriggerEventEnumNet	FirmwareDestinationNames, 691
Mcs::Usb, 63	Mcs::Usb, 53
DigitalStimulatorTriggerSlopeEnumNet	DSPAnalogGroup
Mcs::Usb, 64	Mcs::Usb, 71, 82
DigitalTargetEnumNet	DSPDataGroup
Mcs::Usb, 64	Mcs::Usb, 60, 92
Digout	DSPDigitalGroup
Mcs::Usb, 64	Mcs::Usb, 71, 82
DigOutStim	DummyCommand
Mcs::Usb, 63	CLIH3DeviceNet, 197
DigOutStimulatorStartTrigger	Duration
Mcs::Usb, 64	CStimulusFunctionNet::SidebandData, 711
DigOutStimulatorStopTrigger	,
Mcs::Usb, 64	eCube
Digstream	Mcs::Usb, 78, 82
Mcs::Usb, 64	eCubeHeadstage
DigStreamFromReceiver	Mcs::Usb, 67
Mcs::Usb, 93	ElectricalStimulation
DigStreamToReceiver	HeadStageIDType, 694
Mcs::Usb, 64	ElectrodeDacMuxEnumNet
Dilutor	Mcs::Usb, 64
Mcs::Usb, 79	ElectrodeModeEnumNet
DisableAutoReset	Mcs::Usb, 64
CStg200xDownloadBasicNet, 561	ElectrodeOffset
DisableMultiFileMode	Mcs::Usb, 90
CStg200xDownloadNet, 569	emAutomatic
Joig Look Sommodulion, Ood	Mcs::Usb, 65

emManual	HeadStageIDType, 695
Mcs::Usb, 65	EOFAndCRC
EmptyKey	Mcs::Usb, 60
CMcsUsbNet, 324	Equals
Emu_GetCellCapacity	CMcsUsbListEntryNet, 311
CRoboDacqNet, 460	HeadStageIDType, 695
Emu_GetCellPotential	HeadstageIDTypeObject, 697
CRoboDacqNet, 460	EraseEepromRegisterPreconfig
Emu_GetCellResists	CMcsUsbNet, 324
CRoboDacqNet, 460	EraseFilterParameterPermanent
Emu_GetElectrodeResists	CFilterConfigurationNet, 138
CRoboDacqNet, 460	CFilterConfigurationRegisterNet, 140
Emu_GetNoise	ErasePermanentAdcOffset
CRoboDacqNet, 460	CLIH3DeviceNet, 198
Emu_SetCellCapacity	ErasePermanentDacOffset
CRoboDacqNet, 460	CLIH3DeviceNet, 198
Emu_SetCellPotential	Error_Callback_Aquisition_Stopped
CRoboDacqNet, 460	CMcsUsbDacqNet, 296
Emu_SetCellResists	Error_Callback_Data_lost
CRoboDacqNet, 460	CMcsUsbDacqNet, 296
Emu_SetElectrodeResists	Error_Callback_Frames_Lost
CRoboDacqNet, 460	CMcsUsbDacqNet, 296
Emu_SetNoise	Error_Callback_Packet_Error
CRoboDacqNet, 460	CMcsUsbDacqNet, 296
EnableAnalogOut	Error_Callback_Queue_Full
CSCUFunctionNet, 505	CMcsUsbDacqNet, 297
EnableAutoReset	Error_Callback_RingQueue_Full
CStg200xDownloadBasicNet, 561	CMcsUsbDacqNet, 297
EnableChannelsInGroup	ErrorEvent
CCMOSMea_FunctionNet, 104	CMcsUsbDacqNet, 297
CDacqGroupChannelSelectionTemplateNet< Dac-	ExtendedMultiFileMode
qGroupChannelEnumTemplateNet, Dac-	Mcs::Usb, 84
qGroupChannelEnumTemplate, CDevice-	ExternBCTester
GroupChannelInfoTemplateNet >, 122	Mcs::Usb, 78
EnableChecksum	ExternDTester
CMeaDeviceNet, 355	Mcs::Usb, 78
COctoPotDeviceNet, 398	ExternSTester
EnableDigitalIn	Mcs::Usb, 78
CMeaDeviceNet, 356	FactoryReset
COctoPotDeviceNet, 398	CTcxDeviceNet, 589
EnableExceptions	Falling
CMcsUsbNet, 324	Mcs::Usb, 64
EnableMultiFileMode	FCB
CStg200xDownloadNet, 569	Mcs::Usb, 77
EnableQueue	FCX
CRoboDeviceNet, 475	Mcs::Usb, 77
EnableTimestamp	Feedback
CMeaDeviceNet, 357	Mcs::Usb, 57, 71, 73, 74, 83, 86, 93
COctoPotDeviceNet, 398	FeedbackGetSampleTimerCount
EnableUserTrigger	CMeaFeedbackFunctionNet, 364
CLIH3DeviceNet, 197 Encapsulator	FeedbackGetTriggerTotzeitFactor
Mcs::Usb, 79	CMeaFeedbackFunctionNet, 364
	FeedbackHigh
enCMosMeaChipType Mcs::Usb, 65	Mcs::Usb, 63
EnSTG200x_STATUS	FeedbackLow
Mcs::Usb, 65	Mcs::Usb, 63
Entry	FeedbackReg
Linu y	Mcs::Usb, 59

FeedbackSetAnalogSource	FilterType
CMeaFeedbackFunctionNet, 365	CFilterPropertyNet, 142
FeedbackSetChannelFilter	FilterTypeEnumNet
CMeaFeedbackFunctionNet, 365	Mcs::Usb, 66
FeedbackSetCreateFilter	FindEndpoints
CMeaFeedbackFunctionNet, 365	CGenericDevelopDeviceNet, 158
FeedbackSetDigitalMapping	FindFilter
CMeaFeedbackFunctionNet, 365	CCreateFilterNet, 117
FeedbackSetFeedback	FindFirmwareVersionMagicInBuffer
CMeaFeedbackFunctionNet, 365	CMcsUsbFactoryNet, 302
FeedbackSetFilterOff	FindReference
CMeaFeedbackFunctionNet, 365	CRoboDeviceNet, 475
FeedbackSetFilterParameter	FindReferencel
CMeaFeedbackFunctionNet, 365	CRoboStatorDeviceNet, 494
FeedbackSetFilterParameter32	FindReferencePhase0
CMeaFeedbackFunctionNet, 365	CRoboDeviceNet::RoboMainLowLevelCommands,
FeedbackSetGlobalChannelFilter	703
CMeaFeedbackFunctionNet, 366	FindReferencePhase0XY
FeedbackSetIIRFilterParameter	CRoboStatorDeviceNet::RoboMainStatorLowLevelCommands,
CMeaFeedbackFunctionNet, 366	709
FeedbackSetLogic	FindReferenceXY
CMeaFeedbackFunctionNet, 366	CRoboStatorDeviceNet, 494
FeedbackSetMkFilter	FindReferenceZ
CMeaFeedbackFunctionNet, 366	CRoboStatorDeviceNet, 494
FeedbackSetNumberOfLogics	Finished
CMeaFeedbackFunctionNet, 366	Mcs::Usb, 85
FeedbackSetNumberOfRateCounter	FirePressurePulse
	CPPCFunctionNet, 426
CMeaFeedbackFunctionNet, 366 FeedbackSetNumberOfRateDetectors	
	FirmwareDestinationNames, 690
CMeaFeedbackFunctionNet, 366	Altera, 691
FeedbackSetNumberOfSpikeDetectors	Bootstrap, 691
CMeaFeedbackFunctionNet, 367	BUS1_MCSBUS1, 691
FeedbackSetNumberOfTriggers	BUS1_MCSBUS2, 691
CMeaFeedbackFunctionNet, 367	DSP, 691
FeedbackSetRateCounter	FPGA2, 691
CMeaFeedbackFunctionNet, 367	FPGA3, 691
FeedbackSetRateDetector	FPGA4, 691
CMeaFeedbackFunctionNet, 367	FPGA5, 691
FeedbackSetSpikeDetectorThreshold	FPGA6, 691
CMeaFeedbackFunctionNet, 367	MCSBUS1, 692
FeedbackSetTrigger	MCSBUS10, 692
CMeaFeedbackFunctionNet, 367	MCSBUS11, 692
FeedbackSetTriggerTotzeitFactor	MCSBUS12, 692
CMeaFeedbackFunctionNet, 367	MCSBUS13, 692
FilterActive	MCSBUS2, 692
CFilterPropertyNet, 142	MCSBUS3, 692
FilterAttributeEnumNet	MCSBUS4, 692
Mcs::Usb, 65	MCSBUS5, 692
FilterBand	MCSBUS6, 692
CFilterPropertyNet, 142	MCSBUS7, 692
FilterBandEnumNet	MCSBUS8, 693
Mcs::Usb, 65	MCSBUS9, 693
FilterCalculationDirectionEnumNet	MCU1, 693
Mcs::Usb, 66	PIC, 693
FilterFamily	PIC2, 693
CFilterPropertyNet, 142	PIC3, 693
FilterFamilyEnumNet	PIC4, 693
Mcs::Usb. 66	USB, 693

FluidControlDevice	Mcs::Usb, 56
CSafeISDeviceNet, 501	FPGA3_GOLD
ForceStatusEvent	Mcs::Usb, 55
CStg200xDownloadBasicNet, 561	FPGA4
CStimulusFunctionNet, 577	FirmwareDestinationNames, 691
FPGA10	Mcs::Usb, 55
Mcs::Usb, 55	FPGA4 BASE
FPGA10 BASE	Mcs::Usb, <u>56</u>
Mcs::Usb, 56	FPGA4 GOLD
FPGA10 GOLD	Mcs::Usb, 55
Mcs::Usb, 56	FPGA5
FPGA11	FirmwareDestinationNames, 691
Mcs::Usb, 55	Mcs::Usb, 55
FPGA11 BASE	FPGA5 BASE
-	_
Mcs::Usb, 56	Mcs::Usb, 56
FPGA11_GOLD	FPGA5_GOLD
Mcs::Usb, 56	Mcs::Usb, 55
FPGA12	FPGA6
Mcs::Usb, 55	FirmwareDestinationNames, 691
FPGA12_BASE	Mcs::Usb, 55
Mcs::Usb, 56	FPGA6_BASE
FPGA12_GOLD	Mcs::Usb, 56
Mcs::Usb, 56	FPGA6_GOLD
FPGA13	Mcs::Usb, 55
Mcs::Usb, 55	FPGA7
FPGA13_BASE	Mcs::Usb, 55
Mcs::Usb, 56	FPGA7 BASE
FPGA13 GOLD	Mcs::Usb, <u>56</u>
Mcs::Usb, 56	FPGA7 GOLD
FPGA14	Mcs::Usb, 56
Mcs::Usb, 55	FPGA8
FPGA14 BASE	Mcs::Usb, 55
Mcs::Usb, 56	FPGA8 BASE
FPGA14 GOLD	Mcs::Usb, 56
-	FPGA8_GOLD
Mcs::Usb, 56 FPGA15	
	Mcs::Usb, 56
Mcs::Usb, 55	FPGA9
FPGA15_BASE	Mcs::Usb, 55
Mcs::Usb, 56	FPGA9_BASE
FPGA15_GOLD	Mcs::Usb, 56
Mcs::Usb, 56	FPGA9_GOLD
FPGA16	Mcs::Usb, 56
Mcs::Usb, 55	FPGA_BASE
FPGA16_BASE	Mcs::Usb, 56
Mcs::Usb, 56	FPGA_BOOTSTRAP
FPGA16_GOLD	Mcs::Usb, 56
Mcs::Usb, 56	FPGA_GOLD
FPGA2	Mcs::Usb, 55
FirmwareDestinationNames, 691	FPGA_NORMAL
Mcs::Usb, 55	
FPGA2_BASE	FpgaldEnumNet
Mcs::Usb, 56	Mcs::Usb, 66
FPGA2_GOLD	FrameContextGroup
Mcs::Usb, 55	Mcs::Usb, 61
FPGA3	FromIntPtr
FirmwareDestinationNames, 691	StgStatusNet, 712
	_
Mcs::Usb, 55	FromPtr
FPGA3_BASE	StgStatusNet, 712

FullCharge	Cat A and Cura Danizad Data
FullCharge	GetAccelGyroDesiredRate
Mcs::Usb, 70	CW2100_FunctionNet, 617
FullSpeed	GetAccelGyroEnabled
Mcs::Usb, 70	CW2100_FunctionNet, 617
FunkDongleS	GetAccelRange
Mcs::Usb, 78	CW2100_FunctionNet, 617
FX3MCSDataAddress	GetActiveRunningTableNumber
CMcsUsbFactoryNet, 307	CWarnerValveControllerDeviceNet, 651
FX3MCSDataDeviceIdOffset	GetActiveRunningTableNumberEvent
CMcsUsbFactoryNet, 307	CWarnerValveControllerDeviceNet, 668
FX3MCSDataIFB1ImageOffset	GetAdapterCode
CMcsUsbFactoryNet, 307	CMealmpedanceDeviceNet, 368
FX3MCSDataIFB2ImageOffset	CTEERFunctionNet, 602
CMcsUsbFactoryNet, 307	GetAdapterType
FX3MCSDataVersionOffset	CMcsUsbDacqNet, 269
CMcsUsbFactoryNet, 307	GetAdc
FYIProgram	CFluidControlDeviceNet, 144
•	GetAdcDataFormat
CFYIDeviceNet, 150	
FYITemp	CMcsUsbDacqNet, 270
CFYIDeviceNet, 150	GetADCGain
Gain	CGrapheneFunctionNet, 173, 174
	GetADCInputOffset
CMeaDeviceNet, 360	CCMOSMea_FunctionNet, 105
Gain100	GetADCOffset
Mcs::Usb, 85	CGrapheneFunctionNet, 174
Gain1000	GetAdcOffset
Mcs::Usb, 85	CLIH3DeviceNet, 198
Gain200	COctoPotDeviceNet, 398
Mcs::Usb, 85	GetAdcOffsetU1
Gain500	CTEERFunctionNet, 602
Mcs::Usb, 85	GetAdcOffsetU2
Gate	CTEERFunctionNet, 602
Mcs::Usb, 80	GetADCs
Gated_High_Active	
Mcs::Usb, 79	CMcsBus_SensorNet, 230
Gated_Low_Active	GetADCsLoop
Mcs::Usb, 79	CMcsBus_SensorNet, 231
GE2100	GetAdcZero
	CMcsUsbDacqNet, 270
Mcs::Usb, 78	GetAirpressure
Generic	CRoboDeviceNet, 475
Mcs::Usb, 77	GetAirpressureLimit
Get2AnalogInput	CRoboDeviceNet, 475
CMcsBus_SensorNet, 230	GetAirValve
Get2DigitalInput	CRoboDeviceNet, 475
CMcsBus_SensorNet, 230	GetAllDigout
Get4ADC	CRoboDacqNet, 461
CMcsBus_SensorNet, 230	GetAmplification
Get4ADCAverage	CPgaDeviceNet, 410
CMcsBus_SensorNet, 230	GetAmplitude_nA
Get4ADCCatchampAverageShift	• —
CMcsBus_SensorNet, 230	CTEERFunctionNet, 602
Get4ADCMode	GetAnalogGain
	CMeaDeviceNet, 357
CMcsBus_SensorNet, 230	GetAnalogOutADCRange
Get4DAC	CSCUFunctionNet, 506
CMcsBus_SensorNet, 230	GetAnalogOutChannel
GetAbsMaxCurrentInMicroAmp	CW2100_FunctionNet, 617
CMultiwellOptoStimFunctionNet, 392	GetAnalogOutChannels
GetAccelGyroCurrentRate	CSCUFunctionNet, 506
CW2100_FunctionNet, 617	

GetAnalogOutDACRange	GetAvailableStateListEx
CSCUFunctionNet, 506	CRFFunctionNet, 450
GetAnalogOutFilter	GetAvailableVoltageRangesInMicroVolt
CW2100_FunctionNet, 617	CMcsUsbDacqNet::CHWInfo, 187
GetAnalogRanges	GetAvailableVoltageRangesInMicroVoltAndStringsInMilliVolt
CStg200xBasicNet, 526	CMcsUsbDacqNet::CHWInfo, 187
GetAnalogResolution	GetAxisConfig
CStg200xBasicNet, 526	CRoboDeviceNet::RoboMainLowLevelCommands,
GetAnalogThresholdHigh	703
CWarnerValveControllerDeviceNet, 651	GetAxisLED
GetAnalogThresholdHighEvent	CRoboocyte2DeviceNet, 492
CWarnerValveControllerDeviceNet, 668	GetAxisParametersSignedEeprom
GetAnalogThresholdLow	CMcsBus_AxisParametersNet, 205
CWarnerValveControllerDeviceNet, 652	GetAxisParametersUnsignedEeprom
GetAnalogThresholdLowEvent	CMcsBus_AxisParametersNet, 205
CWarnerValveControllerDeviceNet, 668	GetBaseFrequency
GetAnalogValueUnit	CRFFunctionNet, 451
CMcsUsbDacqNet, 270	GetBaseSamplerate
GetAnalogVoltage	CCMOSMeaDeviceNet, 114
CPPCFunctionNet, 426	GetBath
CPPS_FunctionNet, 435	CCMOSMea_FunctionNet, 105
CWarnerValveControllerDeviceNet, 652	GetBathMode
GetAnalogVoltageEvent	CCMOSMea_FunctionNet, 105
CWarnerValveControllerDeviceNet, 668	GetBatteryState
GetAnalogVoltageRange	CW2100_FunctionNet, 618
CPPCFunctionNet, 426	GetBatteryVoltage
GetAnalogVoltages	CMultiBatteryChargerDeviceNet, 377
CPPS_FunctionNet, 435	GetBiQuad
GetArraySize	
•	CCreateFilterNet, 117
CotAudioChannela	GetBiQuads
GetAudioChannels	CCreateFilterNet, 117
CM9100 FunctionNet, 342	GetBlankingEnable
CW2100_FunctionNet, 617	CStg200xBasicNet, 527
GetAudioOutDacParameter	GetBoardTemp
CLIH3DeviceNet, 198	CTcxDeviceNet, 589
GetAutocalibrationDisabled	GetBoostAlwaysOnMode
CStg200xBasicNet, 526	
	CW2100_StimulatorFunctionNet, 624
GetAvailableBaseSamplerates	GetBoostPreTime
CCMOSMeaDeviceNet, 114	GetBoostPreTime CW2100_StimulatorFunctionNet, 624
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetByteBuffer CGenericDevelopDeviceNet, 158 GetBytesAvailable CSerialPortNet, 520
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519 GetAvailableMemory	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetByteBuffer CGenericDevelopDeviceNet, 158 GetBytesAvailable
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519 GetAvailableMemory CStg200xBasicNet, 527	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetByteBuffer CGenericDevelopDeviceNet, 158 GetBytesAvailable CSerialPortNet, 520
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519 GetAvailableMemory CStg200xBasicNet, 527 CStimulusFunctionNet, 577	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteSusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetBytesAvailable CSerialPortNet, 520 GetBytesPerSample
CCMOSMeaDeviceNet, 114 CGrapheneASICDeviceNet, 170 GetAvailableChambers CWarnerUssingFunctionNet, 633 GetAvailableDeviceList CRFFunctionNet, 450 GetAvailableDeviceListEx CRFFunctionNet, 450 GetAvailableHeadstages CSCUFunctionNet, 506 CW2100_FunctionNet, 617 GetAvailableHeadstagesEvent CSCUFunctionNet, 519 GetAvailableMemory CStg200xBasicNet, 527 CStimulusFunctionNet, 577 GetAvailableSampleRates	GetBoostPreTime CW2100_StimulatorFunctionNet, 624 GetBubbleState CPPS_FunctionNet, 435 GetBubbleStatus CMcsBus_SensorNet, 231 GetBuffer CGenericDevelopDeviceNet, 158 GetBusAddress CMcsBusNet, 244 GetBusAddressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetByteSdressEeprom CMcsBusNet, 244 GetByteBuffer CGenericDevelopDeviceNet, 158 GetBytesAvailable CSerialPortNet, 520 GetBytesPerSample CTEERFunctionNet, 603

OT D : N 500	0144
CTcxDeviceNet, 589	CWarnerUssingFunctionNet, 633
GetCalibrationMax	GetCMOSDataDictionary CCMOSMeaDeviceNet, 115
CTcxDeviceNet, 589	•
GetCalibrationMin	GetCoilCommunication CPositionIIDeviceNet, 413
CTcxDeviceNet, 589 GetCanDisableCalibration	
	GetColorRgb
CStg200xBasicNet, 528 GetCanDisableCurrentswitch	CMultiwellOptoStimFunctionNet, 393 GetColorStr
CStg200xBasicNet, 528	CMultiwellOptoStimFunctionNet, 393
GetCapacityC	GetCommand
CRoboDacqNet, 461	CMcsBusNet, 244, 245
GetCapacityV	CPedoterDeviceNet, 407
CRoboDacqNet, 461	CRoboDacqNet, 461
GetCapacityX	GetComplianceVoltageThreshold
CRoboDacqNet, 461	CWarnerUssingFunctionNet, 633
GetCardinalDacqSamplerate	GetConfiguration
CInterfaceboardFunctionNet, 194	CMcsUsbNet, 325
GetCardinalStgOutputrate	GetConfigurationBit
CInterfaceboardFunctionNet, 194	CRoboDacqNet, 461
GetChannel	GetConfigurationBitAxc
CSw2to64DeviceNet, 585	CRoboDacqNet, 461
GetChannelDataFillSize	GetConfigurationBitBlu_Led
CMcsUsbDacqNet, 270	CRoboDacqNet, 461
GetChannelDatal16	GetConfigurationBitBlu_LedToggleFast
CCMOSMeaDeviceNet, 114	CRoboDacqNet, 461
GetChannelDatal32	GetConfigurationBitBlu_LedToggleSlow
CCMOSMeaDeviceNet, 114	CRoboDacqNet, 461
GetChannelDataUI16	GetConfigurationBitCC_Gen
CCMOSMeaDeviceNet, 114	CRoboDacqNet, 462
GetChannelDataUI32	GetConfigurationBitCV_Gen
CCMOSMeaDeviceNet, 115	CRoboDacqNet, 462
GetChannelLayout	GetConfigurationBitRC_Gen
CMcsUsbDacqNet, 270	CRoboDacqNet, 462
GetChannels	GetConfigurationBitRed_Led
CMultiBatteryChargerDeviceNet, 377	CRoboDacqNet, 462
CSw2to64DeviceNet, 585	GetConfigurationBitRed_LedSaturation
GetChannelsCountOfChamber	CRoboDacqNet, 462
CWarnerUssingFunctionNet, 633	GetConfigurationBitRed_LedToggleFast
GetChannelsInBlock	CRoboDacqNet, 462
CMcsUsbDacqNet, 270	GetConfigurationBitRed_LedToggleSlow
GetChannelState	CRoboDacqNet, 462
CMultiBatteryChargerDeviceNet, 377	GetConfigurationBitRelais
GetChargeCapacity	CRoboDacqNet, 462
CMultiBatteryChargerDeviceNet, 377	GetConfigurationBitRV_Gen
GetChargeCurrent	CRoboDacqNet, 462
CMultiBatteryChargerDeviceNet, 378	GetConfigurationBits
GetChargingMode	CRoboDacqNet, 462
CotChargingPCoefficient	GetConfigurationBitStream
GetChargingPCoefficient CMultiBatteryChargerDeviceNet, 378	CRoboDacqNet, 462 GetConfigurationBitSupply
GetChecksumFromFX3Image	CRoboDacqNet, 463
CMcsUsbFactoryNet, 302	GetConnectedDevice
GetCheckVoltage	CRFFunctionNet, 451
COkuvisionStimulatorDeviceNet, 402	GetControllerParams
GetClampAmpSerialNumber	CTEERFunctionNet, 603
CRoboDacqNet, 461	GetCpp
GetClampMode	CCreateFilterNet, 117
CTEERFunctionNet, 603	GetCrossTalkOffset

CRoboDacqNet, 463	GetDACOffset
GetCrossTalkOptimum	CGrapheneFunctionNet, 175
CRoboDacqNet, 463	COkuvisionStimulatorDeviceNet, 402
GetCur2VolResistance	GetDacOffset
CGrapheneFunctionNet, 175	CDacCalibrationFunctionNet, 119
GetCurrent	CLIH3DeviceNet, 199
CTcxDeviceNet, 589	COctoPotDeviceNet, 398
GetCurrentAirvalve	GetDacPampsPerDigitHighCurrentRange
CRoboDeviceNet, 475	CWarnerUssingFunctionNet, 634
GetCurrentAirvalveLimit	GetDacPampsPerDigitLowCurrentRange
CRoboDeviceNet, 475	CWarnerUssingFunctionNet, 634
GetCurrentCycle	GetDacqRunStatus
CMeaCoatDeviceNet, 349	CLIH3DeviceNet, 199
GetCurrentEditTableNumber	GetDacRange
CWarnerValveControllerDeviceNet, 652	CW2100_FunctionNet, 618
GetCurrentEnable	GetDACResolution
CTEERFunctionNet, 603	CStg200xBasicNet, 530
GetCurrentFactor	CStimulusFunctionNet, 578
COkuvisionStimulatorDeviceNet, 402	CW2100_StimulatorFunctionNet, 624
GetCurrentNumberOfValves	GetDACs
CWarnerValveControllerDeviceNet, 652	CMcsBus_SensorNet, 231
GetCurrentNumberOfValvesEvent	GetDacUseldleValue
CWarnerValveControllerDeviceNet, 669	CLIH3DeviceNet, 199
GetCurrentPosition	GetDacZero
CRoboDeviceNet, 476	CTEERFunctionNet, 604
GetCurrentPositionI	CWarnerUssingFunctionNet, 634
CRoboStatorDeviceNet, 494	GetDataFormat
GetCurrentPositionXY	CMcsUsbDacqNet, 271
CRoboStatorDeviceNet, 494	GetDataMode
GetCurrentPositionZ	CMcsUsbDacqNet, 271
CRoboStatorDeviceNet, 495	GetDDecp
GetCurrentRangeInNanoAmp	CTcxDeviceNet, 590
CStg200xBasicNet, 528	GetDebugData
CStimulusFunctionNet, 577	CPositionIIDeviceNet, 413
CW2100_StimulatorFunctionNet, 624	GetDefaultFilterAttributes
GetCurrentRangeInNanoAmpByIndex	CFilterConfigurationNet, 138
CStg200xBasicNet, 528	GetDestination
GetCurrentRangeListInNanoAmp	CMcsUsbFactoryNet, 302
CStg200xBasicNet, 529	GetDestinationCode
GetCurrentRangeSelectedIndex	DriverVersionNet, 686
CStg200xBasicNet, 529	GetDestinationDisplayLabel
GetCurrentResolutionInNanoAmp	CMcsUsbFactoryNet, 302
CStg200xBasicNet, 529	GetDestinationName
CStimulusFunctionNet, 577	CMcsUsbFactoryNet, 302, 303
CW2100_StimulatorFunctionNet, 624	DriverVersionNet, 687
GetCurrentResolutionInNanoAmpByIndex	GetDestinationSerialNumber
CStg200xBasicNet, 530	CMcsUsbFactoryNet, 303
GetCycle	GetDestinationTargetAddress
CMeaCleanDeviceNet, 345	CMcsUsbFactoryNet, 303
GetCycles	GetDestinationTargetSelectionInfo
CMeaCleanDeviceNet, 345	CMcsUsbFactoryNet, 303
CMeaCoatDeviceNet, 349	GetDetectionThreshold
GetD	CMcsBus_SensorNet, 231
CTcxDeviceNet, 590	GetDetectorValue
GetDacAmplificationFactor	CMcsBus_SensorNet, 231
CStg200xBasicNet, 530	GetDevice
GetDacIdleValue	CTcxDeviceNet, 590
CLIH3DeviceNet, 199	GetDeviceCannotStallOutRequests

CMcsUsbNet, 325	GetDigoutValue
GetDeviceCapableSpeed	CStg200xBasicNet, 531
CMcsUsbNet, 325	GetDIO
•	
GetDeviceEnum	CMcsBus_FYIExtensionNet, 208
CMcsUsbNet, 325	GetDischargeCapacity
GetDeviceGroupChannelInfos	CMultiBatteryChargerDeviceNet, 378
CDacqGroupChannelSelectionTemplateNet< Dac-	GetDischargeCurrent
qGroupChannelEnumTemplateNet, Dac-	CMultiBatteryChargerDeviceNet, 379
qGroupChannelEnumTemplate, CDevice-	GetDischargeCurrentSetPoint
GroupChannelInfoTemplateNet $>$, 122	CMultiBatteryChargerDeviceNet, 379
GetDeviceId	GetDisplayMode
CMcsUsbNet, 325	CWarnerValveControllerDeviceNet, 653
GetDeviceList	GetDisplayModeEvent
CPositionImpDeviceNet, 421	CWarnerValveControllerDeviceNet, 669
GetDeviceNames	GetDisplayText
CRadioControledDevicesNet, 444	CRoboDacqNet, 463
GetDeviceRootHubVendorEnum	GetDMax
CMcsUsbNet, 325	CTcxDeviceNet, 590
GetDeviceRootHubVendorID	GetDMin
CMcsUsbNet, 325	CTcxDeviceNet, 590
GetDeviceRootHubVendorName	GetDownsampleFactor
CMcsUsbNet, 325	CRoboDacqNet, 463
GetDeviceSpeed	GetDSPHighPassByIndex
CMcsUsbNet, 326	CIntanMea_FunctionNet, 190
GetDeviceType	GetDuration
CTcxDeviceNet, 590	CMeaCoatDeviceNet, 350
GetDevname	GetDuty
CTcxDeviceNet, 590	CTcxDeviceNet, 590
GetDigin	GetEEpromPage
CFluidControlDeviceNet, 145	CLIH3DeviceNet, 200
GetDigInState	GetElectrodeDacMux
CLIH3DeviceNet, 200	CStg200xBasicNet, 531
GetDiginValue	GetElectrodeEnable
CStg200xBasicNet, 530	CStg200xBasicNet, 532
GetDigitalData	GetElectrodeMode
CMeaDigitalDataFunctionNet, 362	CStg200xBasicNet, 533
GetDigitalIn	GetEnableAmplifierProtectionSwitch
CPPCFunctionNet, 428	CStg200xBasicNet, 534
CPPS FunctionNet, 435	GetEnabledChannelsInGroup
-	•
GetDigitalOutPortValve	CCMOSMea_FunctionNet, 105
CWarnerValveControllerDeviceNet, 653	CDacqGroupChannelSelectionTemplateNet < Dac
GetDigitalOutPortValveEvent	qGroupChannelEnumTemplateNet, Dao
CWarnerValveControllerDeviceNet, 669	qGroupChannelEnumTemplate, CDevice
GetDigitalPortDirection	GroupChannelInfoTemplateNet >, 122, 123
CWarnerValveControllerDeviceNet, 653	GetEnableHeaterLimit
GetDigitalPortDirectionEvent	CTcxDeviceNet, 591
CWarnerValveControllerDeviceNet, 669	GetEnableThermocouple
GetDigitalSource	CTcxDeviceNet, 591
CMcsUsbDacqNet, 271–273	GetEntry
GetDigitalStimulatorTrigger	CMcsUsbListEntryNet, 312
CW2100_StimulatorFunctionNet, 624	GetEntryCount
GetDigitalStimulatorTriggerSlope	CMcsUsbListEntryNet, 313
CW2100_StimulatorFunctionNet, 625	GetEnumerationSpeed
GetDigout	CMeaDeviceNet, 357
CFluidControlDeviceNet, 145	GetErrorAirpressure
CRoboDacqNet, 463	CRoboDeviceNet, 476
GetDigoutMode	GetErrorCurrentAirvalve
CStg200xBasicNet, 531	CRoboDeviceNet, 476

GetErrorMessage	GetGilsonDevice
CMcsUsbDacqNet, 274	CRoboocyte2DeviceNet, 492
GetErrorText	GetGlobalRepeat
CMcsUsbNet, 326	CDigOutStimulatorFunctionNet, 129
GetErrorVoltage12V	GetGNDI
CRoboDeviceNet, 476	CCMOSMea_FunctionNet, 105
GetErrorVoltage5V	GetGroupADCBits
CRoboDeviceNet, 476	CCMOSMea_FunctionNet, 105
GetErrorVoltageAirvalve	GetGroupChannelBitmaskBySelect
CRoboDeviceNet, 476	CCMOSMea FunctionNet, 105, 106
GetErrorVoltageRs485A	GetGroupChannelBitmaskHS1NCBathCurrent
CRoboDeviceNet, 476	CCMOSMea_FunctionNet, 106
GetErrorVoltageRs485B	GetGroupChannelBitmaskHS1NCCol2Current
CRoboDeviceNet, 476	CCMOSMea_FunctionNet, 106
GetErrorVoltageValves	GetGroupChannelBitmaskHS1NChipTemp
CRoboDeviceNet, 477	CCMOSMea_FunctionNet, 106
GetEventData	GetGroupChannelBitmaskHS1Sidebands
CPositionIIDeviceNet, 414	CCMOSMea FunctionNet, 106, 107
GetExternalElectrodeEnable	GetGroupChannelBitmaskHS1TriggerStatus
CStg200xBasicNet, 535	CCMOSMea_FunctionNet, 107
GetFAAmplification	GetGroupChannelBitmaskIFDigChannels
CStg200xBasicNet, 535	CCMOSMea_FunctionNet, 107
GetFilter	GetGroupChannelBitmaskInterfaceADC
CRoboDacqNet, 463	CCMOSMea_FunctionNet, 107
GetFilterAttributes	GetGroupChannelBitmaskPacketFrameContext
CFilterConfigurationNet, 138	CCMOSMea_FunctionNet, 107, 108
GetFilterCoeffs	GetGroupChannelBitmaskSTG1DACSignal
CRoboDacqNet, 463	CCMOSMea_FunctionNet, 108
GetFilterParametersHeadstage	GetGroupChannelDatal16
-	CMcsUsbDacqNet, 274
CWClassicFunctionNet, 676	·
GetFilterProperties	GetGroupChannelDatal32
CSCUFunctionNet, 506	CMcsUsbDacqNet, 275
CW2100_FunctionNet, 618	GetGroupChannelDataUI16
GetFilterProperty	CMcsUsbDacqNet, 275
CMcsUsbDacqNet, 274	GetGroupChannelDataUI32
CSCUFunctionNet, 507	CMcsUsbDacqNet, 276
CW2100_FunctionNet, 618	GetGroupDCOffset
GetFinalDischargeVoltage	CCMOSMea_FunctionNet, 108
CMultiBatteryChargerDeviceNet, 379	GetGroupID
GetFirmwareVersion	CCMOSMea_FunctionNet, 108
CMcsUsbNet, 326	CDacqGroupChannelSelectionTemplateNet< Dac-
GetFirmwareVersionFromFile	qGroupChannelEnumTemplateNet, Dac-
CMcsUsbFactoryNet, 303	qGroupChannelEnumTemplate, CDevice-
GetFirmwareVersionFromHexFile	GroupChannelInfoTemplateNet >, 123
CMcsUsbFactoryNet, 303	GetGroupNumberOfChannels
GetFPGAFirmwareType	CCMOSMea_FunctionNet, 108, 109
CW2100_FunctionNet, 618	CDacqGroupChannelSelectionTemplateNet< Dac-
GetFrameErrorCounter	qGroupChannelEnumTemplateNet, Dac-
CTEERFunctionNet, 604	qGroupChannelEnumTemplate, CDevice-
GetFrequency	GroupChannelInfoTemplateNet >, 123
CRadioControledDevicesNet, 445	GetGroupResolutionPerDigit
GetFrequencyRange	CCMOSMea_FunctionNet, 109
CPgaDeviceNet, 410	GetGroupSampleSize
GetGain	CCMOSMea_FunctionNet, 109
CMeaDeviceNet, 357	${\tt CDacqGroupChannelSelectionTemplateNet$
CPgaDeviceNet, 411	qGroupChannelEnumTemplateNet, Dac-
GetGate	qGroupChannelEnumTemplate, CDevice-
CCMOSMea FunctionNet, 105	GroupChannelInfoTemplateNet >, 123

GetGroupType CCMOSMea_FunctionNet, 109	CW2100_FunctionNet, 618 CWClassicFunctionNet, 676
CDacqGroupChannelSelectionTemplateNet< Dac-	GetHeadstagePowerStateAtStart
qGroupChannelEnumTemplateNet, Dac-	CSCUFunctionNet, 511
qGroupChannelEnumTemplate, CDevice-	GetHeadstagePresent
GroupChannelInfoTemplateNet >, 123	CMcsUsbNet, 327
GetGroupUnit	GetHeadstageSamplerate
CCMOSMea_FunctionNet, 109	CSCUFunctionNet, 512
GetGyroRange	GetHeadstageSamplingActive
CW2100_FunctionNet, 618	CW2100_FunctionNet, 618
GetHardwareMaxRange	GetHeadstageSerialNumber
CMcsUsbDacqNet, 276	CSCUFunctionNet, 512
GetHardwareMinRange	GetHeaterLimit
CMcsUsbDacqNet, 276	CTcxDeviceNet, 591
GetHardwareRevision	GetHeaterTemp
CMcsUsbNet, 326	CTcxDeviceNet, 591
GetHasChecksum	GetHighCurrentRange
CWClassicFunctionNet, 676	CWarnerUssingFunctionNet, 635
GetHashCode	GetHighpassFilterEnable
HeadstageIDTypeObject, 697	CFilterConfigurationNet, 138
GetHasLowCurrentswitchLatency	GetHWConfig
CStg200xBasicNet, 535	CRoboDeviceNet::RoboMainLowLevelCommands,
GetHasRedLedHeadstage	703
CWClassicFunctionNet, 676	GetHWRevision
GetHasThermocouple	CRoboDeviceNet::RoboMainLowLevelCommands,
CTcxDeviceNet, 591	703
GetHeadstage	GetHWRevisionEeprom
CStg200xBasicNet, 535	CMcsBusNet, 245
GetHeadstageActive	Getl
CMcsUsbNet, 327	CTcxDeviceNet, 591
GetHeadstageAdcBits	GetIC
CSCUFunctionNet, 507	CRoboDacqNet, 463
GetHeadstageAdcRangeInMicroVolt	GetlClamp
CSCUFunctionNet, 507	CRoboDacqNet, 463
GetHeadstageDacBits	GetlCoeff
CSCUFunctionNet, 508	CRobo_FYITemp_FunctionNet, 455
GetHeadstageDacCurrentRangeInMicroAmpere	GetICOffset
CSCUFunctionNet, 508	CRoboDacqNet, 463
GetHeadstageDacCurrentResolutionInNanoAmpere	GetIDecp
CSCUFunctionNet, 508	CTcxDeviceNet, 591
GetHeadstageDacVoltageRangeInMilliVolt	GetIdleModeOffset
CSCUFunctionNet, 509	CWarnerUssingFunctionNet, 635
GetHeadstageDacVoltageResolutionInMicroVolt	GetlGain
CSCUFunctionNet, 509	CRoboDacqNet, 464
GetHeadstageFrameCyclesToCompare	GetIMax
CSCUFunctionNet, 509	CTcxDeviceNet, 592
GetHeadstageGainInPermille	GetIMin
CSCUFunctionNet, 510	CTcxDeviceNet, 592
GetHeadstageID	GetImpedanceResult
CMcsUsbNet, 327	CIntanMea_FunctionNet, 190
CSCUFunctionNet, 510	GetImpedanceTestFrequency
GetHeadstageLinkSpeed	CMealmpedanceDeviceNet, 369
CSCUFunctionNet, 510	GetImpId
GetHeadstageNumberOfAnalogChannels	CPositionImpDeviceNet, 422
CSCUFunctionNet, 511	GetImplantCurrentSetpoint
GetHeadstageNumberOfStimulationChannels	CPositionIIDeviceNet, 414
CSCUFunctionNet, 511	GetImplantResult
GetHeadstageOnOff	CPositionIIDeviceNet, 415

GetImplantState	GetMaxNumberOfHeadstages
CPositionIIDeviceNet, 415	CSCUFunctionNet, 512
GetInMovement	GetMaxNumOfColumns
CRoboDeviceNet, 477	CCMOSMea_FunctionNet, 110
GetIntanRegister	GetMaxP
CIntanMea FunctionNet, 190	CTcxDeviceNet, 592
GetIntBuffer	GetMaxpDecp
CGenericDevelopDeviceNet, 158	CTcxDeviceNet, 592
GetIO	GetMaxpMax
CWarnerValveControllerDeviceTesterFunctionNet,	CTcxDeviceNet, 592
673	GetMaxpMin
GetlOut	CTcxDeviceNet, 593
CTcxDeviceNet, 592	GetMaxPower
GetloVoltage	COkuvisionStimulatorDeviceNet, 402
CInterfaceboard2FunctionNet, 192	CRobo_FYITemp_FunctionNet, 455
GetLastAnswer	GetMaxPressureWaitTime
CGilsonDeviceNet, 169	CRoboDeviceNet::RoboMainLowLevelCommands
GetLastUSBError	703
CMcsUsbNet, 327	GetMaxReadableColumns
GetLatency	CCMOSMeaDeviceNet, 115
CMcsBus SensorNet, 231	GetMaxSamplingFrequency
GetLatencyCounter	CMcsUsbDacqNet, 277
CMcsBus_SensorNet, 231	GetMaxStimulusChannelsPerHeadstage
	_
GetLayoutConfiguration	CSCUFunctionNet, 512
CMEA2100x256FunctionNet, 340	GetMaxVoltage
GetLEDSwitch	CMeaCleanDeviceNet, 345
CMcsBus_ExtensionNet, 207	COkuvisionStimulatorDeviceNet, 402
GetLength	GetMCAcceleration
CRobo_FYIProgram_FunctionNet, 453	CMcsBus_MotorControlNet, 212
GetLiquidResistance	GetMCAccelerationEeprom
CTEERFunctionNet, 604	CMcsBus_MotorControlNet, 212
CWarnerUssingFunctionNet, 635	GetMCAccelerationShortCommand
GetListmodeIndexRange	CMcsBus_MotorControlNet, 212
CStg200xBasicNet, 536	GetMCAxisRevisionEeprom
GetListmodeTriggerSource	CMcsBus_MotorControlNet, 212
CStg200xBasicNet, 536	GetMCBreakCurrent
GetLowCurrentRange	CMcsBus_MotorControlNet, 212
CWarnerUssingFunctionNet, 636	GetMCBreakCurrentEeprom
GetLowerFrequencyByIndex	CMcsBus_MotorControlNet, 213
CIntanMea_FunctionNet, 190	GetMCConfig
GetMajor	CMcsBus_MotorControlNet, 213
DriverVersionNet, 687	GetMCConfigEeprom
GetMaxChunkSize_Byte	CMcsBus_MotorControlNet, 213
CTEERFunctionNet, 604	GetMCCurrent
GetMaxCurrent	CMcsBus_MotorControlNet, 213
CMeaCoatDeviceNet, 350	GetMCCurrentEeprom
GetMaxDurationHighCurrentInMicroSec	CMcsBus_MotorControlNet, 213
CMultiwellOptoStimFunctionNet, 393	GetMCCurrentMode
GetMaxDutyCycleHighCurrent	CMcsBus_MotorControlNet, 213
CMultiwellOptoStimFunctionNet, 393	GetMCCurrentModeEeprom
GetMaxHeaterPowerMultiwell	CMcsBus_MotorControlNet, 213
CTcxDeviceNet, 592	GetMCCurrentModeShortCommand
GetMaxNoPressure	CMcsBus_MotorControlNet, 214
CRoboDeviceNet::RoboMainLowLevelCommands,	GetMCCurrentPosition
703 CatMayNeProceuraWaitTime	CMcsBus_MotorControlNet, 214
GetMaxNoPressureWaitTime	GetMCCurrentShortCommand
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_MotorControlNet, 214
703	GetMCCurrentSpeed

CMcsBus_MotorControlNet, 214	CMcsBus_MotorControlNet, 218
GetMCMaxAcceleration	GetMCStandbyTimeEeprom
CMcsBus_MotorControlNet, 214	CMcsBus MotorControlNet, 218
GetMCMaxAccelerationEeprom	GetMea21UsbPort
CMcsBus_MotorControlNet, 214	CMcsUsbNet, 328
GetMCMaxCurrent	GetMeaLayout
CMcsBus_MotorControlNet, 214	CMcsUsbDacqNet, 277
GetMCMaxCurrentEeprom	GetMemoryUsageDAC
CMcsBus MotorControlNet, 215	CStg200xDownloadBasicNet, 561
GetMCMaxSpeed	GetMemoryUsageSyncout
CMcsBus MotorControlNet, 215	CStg200xDownloadBasicNet, 562
GetMCMaxSpeedEeprom	GetMinimalThreshold
CMcsBus_MotorControlNet, 215	CMcsBus_SensorNet, 231
GetMCMaxTravel	GetMinNoPressureWaitTime
CMcsBus_MotorControlNet, 215	CRoboDeviceNet::RoboMainLowLevelCommands,
GetMCMaxTravelEeprom	703
CMcsBus_MotorControlNet, 215	GetMinor
GetMCMaxTravelShortCommand	DriverVersionNet, 688
CMcsBus_MotorControlNet, 215	GetMinPressure
GetMCMovement	
	CRobo DeviceNet, 477
CMcsBus_MotorControlNet, 215 GetMCNewPosition	CRoboDeviceNet::RoboMainLowLevelCommands, 703
	GetMinPressureWaitTime
CMcsBus_MotorControlNet, 216	
GetMCOutputOnOff	CRoboDeviceNet::RoboMainLowLevelCommands,
CMcsBus_MotorControlNet, 216	704
GetMCPhase	GetMinSamplingFrequencyStepsize
CMcsBus_MotorControlNet, 216	CMcsUsbDacqNet, 277
GetMCPhaseOffset	GetMinVoltage
CMcsBus_MotorControlNet, 216	CMeaCleanDeviceNet, 346
GetMCReference	GetModeSelect
CMcsBus_MotorControlNet, 216	CPulseGeneratorFunctionNet, 442
GetMCReferenceCurrent	GetModuleCurrent
CMcsBus_MotorControlNet, 216	CStg200xDownloadNet, 569
GetMCReferenceCurrentEeprom	GetModuleTemp
CMcsBus_MotorControlNet, 216	CStg200xDownloadNet, 569
GetMCRegulatorGain	GetMovementError
CMcsBus_MotorControlNet, 217	CRoboDeviceNet, 477
GetMCRegulatorGainEeprom	GetMovePump
CMcsBus_MotorControlNet, 217	CMcsBus_SensorNet, 232
GetMcsBus_Extension	GetMultiHeadstageMode
CRoboocyte2DeviceNet, 492	CW2100_FunctionNet, 618
GetMCScalingFactor	GetMultiplexedDataChannelsInBlock
CMcsBus_MotorControlNet, 217	CStimulusFunctionNet, 578
GetMCScalingFactorEeprom	GetNanoVoltsPerKelvin
CMcsBus_MotorControlNet, 217	CMcsBus_TempSensorNet, 237
GetMCSpeed	GetNeurochipMemoryData
CMcsBus_MotorControlNet, 217	CCMOSMea FunctionNet, 110
GetMCSpeedEeprom	GetNeurochipMemorySize
CMcsBus_MotorControlNet, 217	CCMOSMea_FunctionNet, 110
GetMCSpeedShortCommand	GetNIC_MS
CMcsBus_MotorControlNet, 217	CRoboDacqNet, 464
GetMCSpeedUnitEeprom	GetNUC_MS
CMcsBus_MotorControlNet, 218	CRoboDacqNet, 464
GetMCStandbyCurrent	GetNumAmplifications
CMcsBus_MotorControlNet, 218	CPgaDeviceNet, 411
GetMCStandbyCurrentEeprom	GetNumber
CMcsBus_MotorControlNet, 218	CMeaSwitchDeviceNet, 371
GetMCStandbyTime	CSw2to64DeviceNet, 585
activicotationy fillie	OUWELOUTDEVICEINEL, JOJ

GetNumberOfAnalogChannels	GetNUV MS
CStg200xBasicNet, 536	CRoboDacqNet, 464
CStimulusFunctionNet, 578	GetOffsetCurrent
CW2100_StimulatorFunctionNet, 625	CMeaCoatDeviceNet, 350
GetNumberOfAudioChannels	GetOnOff
CMeaAudioFunctionNet, 342	CPositionIIDeviceNet, 415
GetNumberOfAvailableChambers	CTcxDeviceNet, 593
CWarnerUssingFunctionNet, 636	GetOutputCurrent
GetNumberOfAvailableSamples	CMeaCoatDeviceNet, 350
CTEERFunctionNet, 604	GetOutputRate
GetNumberOfChannels	CStg200xBasicNet, 538
CDigOutStimulatorFunctionNet, 130	GetOutputVoltage
GetNumberOfCurrentRangeIndexes	CMeaCleanDeviceNet, 346
CStg200xBasicNet, 537	GetP
GetNumberOfDataBits	CTcxDeviceNet, 593
CMcsUsbDacqNet, 277	GetParameter
GetNumberOfDevices	CRoboDeviceNet::RoboMainLowLevelCommands,
CMcsUsbListNet, 316	704
GetNumberOfHardwareSlotsForChambers	GetPattern
CWarnerUssingFunctionNet, 636	CMeaSwitchDeviceNet, 371
GetNumberOfHWADCChannels	GetPatternBool
CMcsUsbDacqNet::CHWInfo, 188	CMeaSwitchDeviceNet, 371
GetNumberOfHWDACPaths	GetPauseDuration
CStg200xBasicNet, 537	CMeaCoatDeviceNet, 350
GetNumberOfHWDigitalChannels	GetPCoeff
CMcsUsbDacqNet::CHWInfo, 188	CRobo_FYITemp_FunctionNet, 455
GetNumberOfStimulationElectrodes	GetPDecp
CStg200xBasicNet, 537	CTcxDeviceNet, 593
GetNumberOfStimulationSourcesPerElectrode	GetPeriod
CStg200xBasicNet, 537	CPulseGeneratorFunctionNet, 442
GetNumberOfSupportedGroups	GetPeriod_us
CCMOSMea_FunctionNet, 110	CTEERFunctionNet, 605
CDacqGroupChannelSelectionTemplateNet< Dac-	GetPermanentCurrentInMicroAmp
qGroupChannelEnumTemplateNet, Dac-	CMultiwellOptoStimFunctionNet, 394
qGroupChannelEnumTemplate, CDevice-	GetPGain
GroupChannelInfoTemplateNet $>$, 124	CRoboDacqNet, 464
GetNumberOfSyncoutChannels	GetPhases
CStg200xBasicNet, 537	CRoboDeviceNet::RoboMainLowLevelCommands,
CW2100_StimulatorFunctionNet, 625	704
GetNumberOfTriggerInputs	GetPicFirmwareType
CStg200xBasicNet, 537	CW2100_FunctionNet, 618
CW2100_StimulatorFunctionNet, 625	GetPiezoState
GetNumberOfVoltageRangeIndexes	CMcsBus_SensorNet, 232
CStg200xBasicNet, 538	GetPlateClampLockState
GetNumConfigurations	CMultiwellDeviceNet, 386
CMcsUsbNet, 328	GetPlateClampState
GetNumControlChannels	CMultiwellDeviceNet, 386
CTcxDeviceNet, 593	GetPlateClampStateByHeadstage
GetNumDestinations	CMultiwellCallbackFunctionNet, 383
CMcsUsbFactoryNet, 304	GetPlateClampStateByHeadstageEvent
GetNumDevices	CMultiwellCallbackFunctionNet, 383
CTcxDeviceNet, 593	GetPlateMux
GetNumEntries	CMultiwellDeviceNet, 386, 387
DriverVersionNet, 688	GetPlateType
GetNumFrequencyRanges	CMultiwellDeviceNet, 387
CPgaDeviceNet, 411	GetPMax
GetNumMeasureChannels	CTcxDeviceNet, 593
CTcxDeviceNet, 593	GetPMin

CTcxDeviceNet, 593	CSCUFunctionNet, 513
GetPoti	GetReferenceTemperature
CMcsUsbDacqNet, 277	CFluidControlDeviceNet, 145
GetPOut	GetRegulationTimeouts
CTcxDeviceNet, 593	CMcsBus SensorNet, 232
GetPowerMuxPlate	GetRegulatorFactor
CMultiwellDeviceNet, 387	CMcsBus_SensorNet, 233
GetPowerStrength	GetRegulatorOnOff
CPositionIIDeviceNet, 416	CMcsBus_SensorNet, 233
GetPressure	CRobo_FYITemp_FunctionNet, 455
CMcsBus_SensorNet, 232	GetRegulatorStatus
GetPressureOffset	CMcsBus_SensorNet, 233
CMcsBus_SensorNet, 232	GetRepeats
GetPressureRange	CProgramPressureCurveNet, 440
CPPCFunctionNet, 428	GetRes1
GetPulseform	CTcxDeviceNet, 594
COkuvisionStimulatorDeviceNet, 402	GetRes2
GetPulseLength	CTcxDeviceNet, 594
CPulseGeneratorFunctionNet, 442	GetResetFilter
GetPumpCouple	CWClassicFunctionNet, 676
CPPS_FunctionNet, 435	GetResistanceC
GetPumpEnableSpeedRatio	CRoboDacqNet, 464
CPPS FunctionNet, 435	GetResistanceV
GetPumpFastOnOff	CRoboDacqNet, 464
CPPS_FunctionNet, 435	GetResolutionPerDigit
GetPumpFastSpeed	CMcsUsbDacqNet, 278
CPPS_FunctionNet, 435	GetResS
GetPumpFunctionSpeeds	CTcxDeviceNet, 594
CPPS_FunctionNet, 436	GetResult
GetPumpManualOnOff	CMealmpedanceDeviceNet, 369
CPPS_FunctionNet, 436	GetResX
GetPumpMaxSpeed	CTcxDeviceNet, 594
CPPS FunctionNet, 436	GetRFConnectionStatus
GetPumpModeType	CWClassicFunctionNet, 676
CPPCFunctionNet, 428	GetRFFrequency
CPPS_FunctionNet, 436	CPositionImpDeviceNet, 422
GetPumpSpeed	GetRFFrequencyHeadstage
CRoboFluidDeviceNet, 488	CWClassicFunctionNet, 676
GetPumpSpeedRatio	GetRFFrequencyReceiver
• •	·
CPPS_FunctionNet, 436	CWClassicFunctionNet, 676
GetPumpSpeedUnit	GetRFPower
CPPCFunctionNet, 429	CWClassicFunctionNet, 677
CPPS_FunctionNet, 436	GetRoboDacq
GetPWM	CRoboocyte2DeviceNet, 492
CFluidControlDeviceNet, 145	GetRoboFluidDevice
GetPwrOut	CEncapsulatorDeviceNet, 133
CTcxDeviceNet, 594	CRoboocyte2DeviceNet, 492
GetPwrSet	GetRotaryPositionCode
CTcxDeviceNet, 594	CTEERFunctionNet, 605
GetRatedCapacity	GetRotatePump
CMultiBatteryChargerDeviceNet, 380	CMcsBus_SensorNet, 233
GetReady	GetROut
CMealmpedanceDeviceNet, 369	CTcxDeviceNet, 594
GetRecordingNumber	GetRTC
CRoboDacqNet, 464	COkuvisionStimulatorDeviceNet, 403
GetReferenceElectrodeMode	CPositionIIDeviceNet, 416
CSCUFunctionNet, 513	GetSampleBufferChunk
GetReferenceElectrodeSwitchState	CTEERFunctionNet, 605

CatCamplaIntanual	CatCatnaintMin
GetSampleInterval	GetSetpointMin
CLIH3DeviceNet, 200	CTcxDeviceNet, 595
GetSamplePeriode	GetShortBuffer
CMcsBus_SensorNet, 233	CGenericDevelopDeviceNet, 159
GetSampleRate	GetSimulation
CTEERFunctionNet, 605	CRoboDacqNet, 464
GetSamplerate	GetSingleHeater
CMcsUsbDacqNet, 278	CMcsBus_FYIExtensionNet, 208
GetSampleVoltageBuffer_uV	GetSingleValve
CTEERFunctionNet, 606	CFluidControlDeviceNet, 146
GetScaleFactorU1	CRoboFluidDeviceNet, 488
CTEERFunctionNet, 606	GetSlope
GetScaleFactorU2	CMeaCleanDeviceNet, 346
CTEERFunctionNet, 606	
	CMeaCoatDeviceNet, 351
GetScanHeadstagesResult	GetSoftwareKey
CWClassicFunctionNet, 677	CMcsUsbNet, 328
GetScreen	GetSoftwareKeyString
CRoboDacqNet, 464	CMcsUsbNet, 328
GetSearchReferenceFastAccel	GetSollPressure
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_SensorNet, 233
704	GetSollTemp
GetSearchReferenceFastSpeed	CRobo_FYITemp_FunctionNet, 455
CRoboDeviceNet::RoboMainLowLevelCommands,	GetSourceBulk
704	CCMOSMea_FunctionNet, 110
GetSearchReferenceFineAccel	GetSourceDrain
CRoboDeviceNet::RoboMainLowLevelCommands,	CCMOSMea_FunctionNet, 110
704	GetSourceGate
GetSearchReferenceFineSpeed	CCMOSMea_FunctionNet, 110
CRoboDeviceNet::RoboMainLowLevelCommands,	GetStartTriggerSlope
704	CDigOutStimulatorFunctionNet, 130
GetSearchReferenceMethod	GetState
CRoboDeviceNet::RoboMainLowLevelCommands,	CRFFunctionNet, 451
705	CRobo_FYIProgram_FunctionNet, 454
GetSearchReferenceMoveOut	GetStateDebugData
CRoboDeviceNet::RoboMainLowLevelCommands,	CPositionIIDeviceNet, 416
705	GetStateEventData
GetSearchReferenceOffsetPos	CPositionIIDeviceNet, 417
CRoboDeviceNet::RoboMainLowLevelCommands,	GetStatus
705	CMcsUsbNet, 328
GetSelectedChannels	DriverVersionNet, 689
CW2100_FunctionNet, 619	GetStatusOfLastCommand
GetSelectedHeadstage	CMcsUsbNet, 328
CWClassicFunctionNet, 677	GetStgProgramInfo
	0 0
GetSelectedHeadstageState	CStg200xBasicNet, 538
CW2100_FunctionNet, 619	GetStgVersionInfo
GetSensorType	CStg200xBasicNet, 539
CTcxDeviceNet, 595	GetStimulationPatternMemory
GetSerialNumber	CW2100_StimulatorFunctionNet, 625
CMcsUsbNet, 328	GetStimulatorStatus
DriverVersionNet, 688	COkuvisionStimulatorDeviceNet, 403
GetSerialNumberHeadstage	GetStimulusParametersCache
CWClassicFunctionNet, 677	CW2100_FunctionNet, 619
GetSetpoint	GetStimulusParametersFromSelectedHS
CTcxDeviceNet, 595	CW2100 FunctionNet, 619
GetSetpointDecp	GetStimulusSites
CTcxDeviceNet, 595	CCMOSMea_FunctionNet, 110
	GetStiumlusParameters
GetSetpointMax	
CTcxDeviceNet, 595	CW2100 FunctionNet, 619

GetStopTriggerSlope	CStg200xBasicNet, 540
CDigOutStimulatorFunctionNet, 130	CStimulusFunctionNet, 578
GetSubChannel	GetTotalNumberOfDigitalPorts
CMcsBus_MotorControlNet, 218	CWarnerValveControllerDeviceNet, 654
GetSupplyVoltage	GetTotalNumberOfTables
CPPCFunctionNet, 429	CWarnerValveControllerDeviceNet, 654
CPPS_FunctionNet, 436	GetTotalNumberOfValves
GetSweepCount	CWarnerValveControllerDeviceNet, 654
CStg200xDownloadBasicNet, 562	GetTotalTableSize
GetSync	CWarnerValveControllerDeviceNet, 655
CWarnerValveControllerDeviceTesterFunctionNet,	GetTouchPadEnable
673	CMultiwellDeviceNet, 388
GetSyncoutMap	
CStg200xBasicNet, 539	GetTrigger CStg200xDownloadBasicNet, 562
_	_
GetSyncState	GetTriggerSource
CMcsBus_SensorNet, 233	CStg200xBasicNet, 540
GetTableName	GetU1Offset
CWarnerValveControllerDeviceNet, 653	CWarnerUssingFunctionNet, 636
GetTableNamebyIndex	GetU1Reference
CWarnerValveControllerDeviceNet, 654	CWarnerUssingFunctionNet, 638
GetTableNamebyIndexEvent	GetU2Offset
CWarnerValveControllerDeviceNet, 669	CWarnerUssingFunctionNet, 638
GetTablepointer	GetU2Reference
CRetinaLedDeviceNet, 447	CWarnerUssingFunctionNet, 638
GetTemperatur	GetUByteBuffer
CMcsBus_TempSensorNet, 237	CGenericDevelopDeviceNet, 160
GetTestMode	GetUC
CRFFunctionNet, 451	CRoboDacqNet, 464
GetThermocoupleCalibration	GetUClamp
CFluidControlDeviceNet, 146	CRoboDacqNet, 465
CTcxDeviceNet, 595	GetUCOffset
GetThermocoupleNanovoltPerKelvin	CRoboDacqNet, 465
CFluidControlDeviceNet, 146	GetUintA
CTcxDeviceNet, 595	CFilterCoefficientsNet, 136
GetThermocoupleReferenceTemp	GetUintB
CTcxDeviceNet, 596	CFilterCoefficientsNet, 136
GetThermocoupleTemp	GetUIntBuffer
CTcxDeviceNet, 596	CGenericDevelopDeviceNet, 160
GetThermocoupleTempAbs	GetUnit
CTcxDeviceNet, 596	CTcxDeviceNet, 596
GetThermocoupleTemperature	GetUnitDescription
CFluidControlDeviceNet, 147	CWarnerUssingFunctionNet, 639
GetThermoOffset	GetUnitExponent
CMcsBus_TempSensorNet, 237	CWarnerUssingFunctionNet, 639
GetThermoTemp	GetUnitName
CMcsBus_TempSensorNet, 238	CWarnerUssingFunctionNet, 639
GetThermoVoltage	GetUnitsPerDigit
CMcsBus_TempSensorNet, 238	CWarnerUssingFunctionNet, 640
GetTimeInPause	GetUOut
CMeaCoatDeviceNet, 351	CTcxDeviceNet, 596
GetTimeInPlateau	GetUpdateDisplay
CMeaCoatDeviceNet, 351	CRoboDacqNet, 465
GetTimeResolutionInNanoSeconds	GetUpperFrequencyByIndex
CStg200xBasicNet, 539	CIntanMea_FunctionNet, 190
CW2100_StimulatorFunctionNet, 625	GetUptimeSeconds
GetTimeSlot	CTEERFunctionNet, 606
CW2100_StimulatorFunctionNet, 625	CWarnerUssingFunctionNet, 640
GetTotalMemory	GetUSBDeviceIDFromFX3Image

CMcsUsbFactoryNet, 304	CWarnerValveControllerDeviceNet, 670
GetUsbListEntries	GetValveManualGroup
CMcsUsbListNet, 316	CWarnerValveControllerDeviceNet, 656
GetUsbListEntry	GetValveManualGroupEvent
CMcsUsbListNet, 316	CWarnerValveControllerDeviceNet, 670
CMcsUsbNet, 329	GetValveManualState
GetUseBubble	CWarnerValveControllerDeviceNet, 657
CPPS_FunctionNet, 436	GetValveManualStateEvent
GetUsercodeFromBitFile	CWarnerValveControllerDeviceNet, 670
CMcsUsbFactoryNet, 304	GetValveMode
GetUsercodeFromFlash	CWarnerValveControllerDeviceNet, 657
CMcsUsbFactoryNet, 304	GetValveModeEvent
GetUserDefinedName	CWarnerValveControllerDeviceNet, 670
CW2100_FunctionNet, 619	GetValves
GetUserDefinedNameCache	CMcsBus_FYIExtensionNet, 208
CW2100_FunctionNet, 619	GetValvesActiveMap
GetUserDefinedNameFromSelectedHS	CWarnerValveControllerDeviceNet, 657
CW2100_FunctionNet, 620	GetValvesManualStateMap
GetUserParameter	CWarnerValveControllerDeviceNet, 657
CRoboDeviceNet::RoboMainLowLevelCommands,	GetValveTableEntry
705	CWarnerValveControllerDeviceNet, 658
GetUShortBuffer	GetVDD3I
CGenericDevelopDeviceNet, 161	CCMOSMea_FunctionNet, 111
GetUV	GetVDDI
CRoboDacqNet, 465	CCMOSMea_FunctionNet, 111
GetUVOffset	GetVdsVgs
CRoboDacqNet, 465	CGrapheneFunctionNet, 176
GetValue	GetVdVs
CTcxDeviceNet, 596	CGrapheneFunctionNet, 176, 177
GetValueHires	GetVdVsDAC
CTcxDeviceNet, 597	CGrapheneFunctionNet, 177
GetValve	GetVersion
CFluidControlDeviceNet, 147	CMcsUsbNet, 329
CRoboFluidDeviceNet, 488	GetVersionInt
GetValve1	DriverVersionNet, 689
CRobo_FYIProgram_FunctionNet, 454	GetVersionString
GetValve2	DriverVersionNet, 690
CRobo_FYIProgram_FunctionNet, 454	GetVMMaxNegativeCurrent
GetValveActive	CMcsBus_VoltageModeNet, 240
CPPCFunctionNet, 429	GetVMMaxNegativeCurrentEeprom
CWarnerValveControllerDeviceNet, 655	CMcsBus_VoltageModeNet, 240
GetValveActiveEvent	GetVMMaxNegativeVoltage
CWarnerValveControllerDeviceNet, 669	CMcsBus_VoltageModeNet, 240
GetValveBoardRevision	GetVMMaxNegativeVoltageEeprom
CWarnerValveControllerDeviceNet, 655	CMcsBus_VoltageModeNet, 240
GetValveBoardRevisionEvent	GetVMMaxPositiveCurrent
CWarnerValveControllerDeviceNet, 669	CMcsBus_VoltageModeNet, 240
GetValveBoardRevisionString	GetVMMaxPositiveCurrentEeprom
CWarnerValveControllerDeviceNet, 655	CMcsBus_VoltageModeNet, 240
GetValveCurrent	— · · ·
	GetVMMaxPositiveVoltage
CWarnerValveControllerDeviceNet, 655	CMcsBus_VoltageModeNet, 240
GetValveDigitalInPort	GetVMMaxPositiveVoltageEeprom
CWarnerValveControllerDeviceNet, 656	CMcsBus_VoltageModeNet, 241
GetValveDigitalInPortEvent	GetVMOutputOnOff
CWarnerValveControllerDeviceNet, 670	CMcsBus_VoltageModeNet, 241
GetValveLedOn	GetVMVoltage
CWarnerValveControllerDeviceNet, 656	CMcsBus_VoltageModeNet, 241
GetValveLedOnEvent	GetVolatileClampOffset

CMultiwellDeviceNet, 388	GetVoltageValvesLimit
GetVoltage	CRoboDeviceNet, 478
COkuvisionStimulatorDeviceNet, 403	GetVolti
GetVoltage12V	CTcxDeviceNet, 597
CRoboDeviceNet, 477	GetWaveform
GetVoltage12VLimit	CTEERFunctionNet, 606
CRoboDeviceNet, 477	GetWaveLengthInNanometer
GetVoltage5V	CMultiwellOptoStimFunctionNet, 394
CRoboDeviceNet, 477	GetWorkingFrequency
GetVoltage5VLimit	CRFFunctionNet, 452
CRoboDeviceNet, 477	GetWPADebugMode
GetVoltageAirvalve	CWClassicFunctionNet, 677
CRoboDeviceNet, 477	GetWPAType
GetVoltageAirvalveLimit	CWClassicFunctionNet, 677
CRoboDeviceNet, 477	GetXGain
GetVoltageClampControllerParam_D	CRoboDacqNet, 465
CWarnerUssingFunctionNet, 640	GetXilinxFlashOffset
GetVoltageClampControllerParam_I	CMcsUsbFactoryNet, 304
CWarnerUssingFunctionNet, 641	GetXilinxFlashReadCommand
GetVoltageClampControllerParam_P	CMcsUsbFactoryNet, 304
CWarnerUssingFunctionNet, 641	GND SWITCH BIT
GetVoltageRange	CW2100_StimulatorFunctionNet, 628
CGrapheneFunctionNet, 177, 178	Graphene_ASIC
GetVoltageRangeIndex	Mcs::Usb, 78
CMcsUsbDacqNet, 278	GrapheneASIC
GetVoltageRangeInMicroVolt	Mcs::Usb, 82
CMcsUsbDacqNet, 278	GrapheneASICHeadstage
CStg200xBasicNet, 540	Mcs::Usb, 67
CStimulusFunctionNet, 578	GrapheneProjectTestDevice
CW2100_StimulatorFunctionNet, 625	Mcs::Usb, 79
GetVoltageRangeInMicroVoltByIndex	Ground
CStg200xBasicNet, 540	Mcs::Usb, 57, 64
GetVoltageRangeInMilliVolt	GroupID
CMcsUsbDacqNet, 278	CDeviceGroupChannelInfoTemplateNet< Dacq-
GetVoltageRangeListInMicroVolt	GroupChannelEnumTemplateNet >, 127
CStg200xBasicNet, 541	GroupType
GetVoltageRangeSelectedIndex	CDeviceGroupChannelInfoTemplateNet< Dacq-
CStg200xBasicNet, 541	GroupChannelEnumTemplateNet >, 127
GetVoltageReached	GyroOnly
CGrapheneFunctionNet, 178	Mcs::Usb, 91
GetVoltageResolution	,-
CGrapheneFunctionNet, 178, 179	Hardware
GetVoltageResolutionInMicroVolt	Mcs::Usb, 66
CStg200xBasicNet, 541	HasAnalogOut
CStimulusFunctionNet, 579	CSCUFunctionNet, 513
CW2100 StimulatorFunctionNet, 626	HasDigitalPort
GetVoltageResolutionInMicroVoltByIndex	CStg200xBasicNet, 542
CStg200xBasicNet, 541	HasFilterAmp
GetVoltageRs485A	CStg200xBasicNet, 542
CRoboDeviceNet, 478	HasGalvanicIsolation
GetVoltageRs485ALimit	CSCUFunctionNet, 513
CRoboDeviceNet, 478	HasHSPowerSwitch
GetVoltageRs485B	CSCUFunctionNet, 514
CRoboDeviceNet, 478	HasIMU
GetVoltageRs485BLimit	HeadStageIDType, 695
CRoboDeviceNet, 478	HasOptoCurrentMessurement
GetVoltageValves	HeadStageIDType, 695
CRoboDeviceNet, 478	HasRadioControl
3.1000=01100.100	CRadioControledDevicesNet, 445

HasRef	State, 698
CRoboDeviceNet::RoboMainLowLevelCommands,	HeadstageType
706	HeadStageIDType, 695
HasRefl	HeadstageTypeEnum
CRoboStatorDeviceNet, 495	HeadStageIDType, 694
HasRefXY	HEKA_LIH3_DEVICE
CRoboStatorDeviceNet, 495	Mcs::Usb, 62
HasRefZ	HekaEPC10Double
CRoboStatorDeviceNet, 495	Mcs::Usb, 77
HasSoftwareKey	HekaEPC10Quadro
CMcsUsbNet, 329	Mcs::Usb, 77
Headstage1NCBathCurrentGroup	HekaEPC10Single
Mcs::Usb, 60	Mcs::Usb, 77
Headstage1NCCol2CurrentGroup	HekaEPC10Triple
Mcs::Usb, 60	Mcs::Usb, 77
Headstage1NChipTempGroup	HekaEPCLite
Mcs::Usb, 61	Mcs::Usb, 77
HeadstageElectrodeGroup	HekalTEV100
Mcs::Usb, 60	Mcs::Usb, 77
HeadstageldEnumNet	HekaLIH30
Mcs::Usb, 67	Mcs::Usb, 77
HeadStageIDType, 693	HekaLIH406
CompareTo, 695	Mcs::Usb, 77
ElectricalStimulation, 694	HekaLIH816
Entry, 695	Mcs::Usb, 77
Equals, 695	HekaPG610
HasIMU, 695	Mcs::Usb, 77
HasOptoCurrentMessurement, 695	HekaPG611
HeadStageIDType, 694	Mcs::Usb, 77 HekaPG612
HeadstageType, 695	
HeadstageTypeEnum, 694	Mcs::Usb, 77 HekaPG618
ID, 695	
MeasuringOnly, 694 NumberOfAnalogChannels, 695	Mcs::Usb, 77 HekaPG690
NumberOfAttalogChamels, 695 NumberOfStimulationChannels, 696	Mcs::Usb, 77
OpticalStimulation, 694	HiClamp
SN, 696	Mcs::Usb, 79
StimulusParameters, 696	HiClamp4Uart
ToString, 695	Mcs::Usb, 79
Type, 696	Highpass
TypeValue, 696	Mcs::Usb, 66
Unknown, 694	HighSpeed
UserDefinedName, 696	Mcs::Usb, 70
Valid, 696	HLA
W16lsW14, 696	Mcs::Usb, 77
HeadstageIDTypeObject, 696	HLADacq
AdditionalText, 697	CHLADeviceNet, 186
_IdType, 697	Hs1Digital
AdditionalText, 698	Mcs::Usb, 63
Equals, 697	HS1DigitalData1
GetHashCode, 697	Mcs::Usb, 86
HeadstageIDTypeObject, 697	HS1ElectrodeGroup
ldType, 698	Mcs::Usb, 71
ToString, 697	HS1Sideband1
HeadStageIDTypeState, 698	Mcs::Usb, 57, 72-74
ControlState, 698	HS1Sideband10
DataState, 698	Mcs::Usb, 72
ldType, 698	HS1Sideband11

M 11 1 70	M 111 57 74 70 74
Mcs::Usb, 72	Mcs::Usb, 57, 71, 73, 74
HS1Sideband12	HS1Trigger3Status
Mcs::Usb, 72 HS1Sideband13	Mcs::Usb, 57, 71, 73
	HS1Trigger4Status
Mcs::Usb, 72	Mcs::Usb, 71, 73
HS1Sideband14	HS1Trigger5Status
Mcs::Usb, 72 HS1Sideband15	Mcs::Usb, 71, 73
Mcs::Usb, 72	HS1Trigger6Status Mcs::Usb, 71, 73
HS1Sideband16	HS1Trigger7Status
Mcs::Usb, 72	Mcs::Usb, 71
HS1Sideband17	HS1Trigger8Status
Mcs::Usb, 72	Mcs::Usb, 72
HS1Sideband18	HS1Trigger9Status
Mcs::Usb, 72	Mcs::Usb, 72
HS1Sideband2	Hs2Digital
Mcs::Usb, 57, 72–74	Mcs::Usb, 63
HS1Sideband3	HS2DigitalData1
Mcs::Usb, 57, 72, 73	Mcs::Usb, 86
HS1Sideband4	HS2ElectrodeGroup
Mcs::Usb, 72, 73	Mcs::Usb, 71
HS1Sideband5	HS2Sideband1
Mcs::Usb, 72, 73	Mcs::Usb, 72, 74, 75
HS1Sideband6	HS2Sideband10
Mcs::Usb, 72, 73	Mcs::Usb, 73
HS1Sideband7	HS2Sideband11
Mcs::Usb, 72	Mcs::Usb, 73
HS1Sideband8	HS2Sideband12
Mcs::Usb, 72	Mcs::Usb, 73
HS1Sideband9	HS2Sideband13
Mcs::Usb, 72	Mcs::Usb, 73
Hs1SidebandHigh	HS2Sideband14
Mcs::Usb, 63	Mcs::Usb, 73
Hs1SidebandLow	HS2Sideband15
Mcs::Usb, 63	Mcs::Usb, 73
Hs1Trigger	HS2Sideband16
Mcs::Usb, 63	Mcs::Usb, 73
HS1Trigger10Status	HS2Sideband17
Mcs::Usb, 72	Mcs::Usb, 73
HS1Trigger11Status	HS2Sideband18
Mcs::Usb, 72	Mcs::Usb, 73
HS1Trigger12Status	HS2Sideband2
Mcs::Usb, 72	Mcs::Usb, 73–75
HS1Trigger13Status	HS2Sideband3
Mcs::Usb, 72	Mcs::Usb, 73, 74
HS1Trigger14Status	HS2Sideband4
Mcs::Usb, 72	Mcs::Usb, 73, 74
HS1Trigger15Status	HS2Sideband5
Mcs::Usb, 72	Mcs::Usb, 73, 74
HS1Trigger16Status	HS2Sideband6
Mcs::Usb, 72	Mcs::Usb, 73, 74
HS1Trigger17Status	HS2Sideband7
Mcs::Usb, 72	Mcs::Usb, 73
HS1Trigger18Status	HS2Sideband8
Mcs::Usb, 72	Mcs::Usb, 73
HS1Trigger1Status	HS2Sideband9
Mcs::Usb, 57, 71, 73, 74 HS1Trigger2Status	Mcs::Usb, 73 Hs2SidebandHigh
1101 mggc120tatus	1 1320 luebanui ligit

Mcs::Usb, 63	IFB30GoldenInterfaceboard
Hs2SidebandLow	Mcs::Usb, 67
Mcs::Usb, 63	IFChannel1
Hs2Trigger	Mcs::Usb, 59
Mcs::Usb, 63	IFChannel2
HS2Trigger10Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel3
HS2Trigger11Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel4
HS2Trigger12Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel5
HS2Trigger13Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel6
HS2Trigger14Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel7
HS2Trigger15Status	Mcs::Usb, 59
Mcs::Usb, 72	IFChannel8
•	
HS2Trigger16Status	Mcs::Usb, 59
Mcs::Usb, 72	IFDigChannelsGroup
HS2Trigger17Status	Mcs::Usb, 61, 71, 82, 93
Mcs::Usb, 72	Ignore
HS2Trigger18Status	Mcs::Usb, 80
Mcs::Usb, 72	IM16KRC
HS2Trigger1Status	Mcs::Usb, 79
Mcs::Usb, 72–74	IM16S16KRA
HS2Trigger2Status	Mcs::Usb, 79
Mcs::Usb, 72–74	IM16S8KRA
HS2Trigger3Status	Mcs::Usb, 79
Mcs::Usb, 72, 73	IM64KRB
HS2Trigger4Status	Mcs::Usb, 79
Mcs::Usb, 72, 74	IM64KRC
HS2Trigger5Status	Mcs::Usb, 79
Mcs::Usb, 72, 74	Input
HS2Trigger6Status	Mcs::Usb, 76
Mcs::Usb, 72, 74	Intel
HS2Trigger7Status	Mcs::Usb, 90
Mcs::Usb, 72	InterfaceADCGroup
HS2Trigger8Status	Mcs::Usb, 60, 71, 82, 91
Mcs::Usb, 72	InterfaceBoard2
HS2Trigger9Status	Mcs::Usb, 67
Mcs::Usb, 72	IntToDouble
HWInfo	Mcs::Usb, 66
CMcsUsbDacqNet, 278	InvitroSignalCollectorUnit
HwVersion	Mcs::Usb, 68
CMcsUsbListEntryNet, 314	InvivoSignalCollectorUnit
555552.5t.2t.; j. t6.; 6 t t	Mcs::Usb, 68
ID	IoVoltageEnumNet
HeadStageIDType, 695	Mcs::Usb, 68
Idle	IS32KRA
Mcs::Usb, 85	Mcs::Usb, 79
IdProduct	IsAnalogOutEnabled
DeviceIdNet, 681	CSCUFunctionNet, 514
IdType	
HeadstageIDTypeObject, 698	IsAutomaticAnalogOut
HeadStageIDTypeState, 698	CSCUFunctionNet, 514
IdVendor	IsBusy
DeviceIdNet, 681	CPPCFunctionNet, 429
IFB2GoldenInterfaceboard	IsChamberAvailable
Mcs::Usb, 67	CWarnerUssingFunctionNet, 641
19103030, 01	

IsChipPowered	IsValveOpenEvent
CCMOSMea_FunctionNet, 111	CWarnerValveControllerDeviceNet, 671
IsConnected	IsValveOpenInAnalogMode
CMcsUsbNet, 329	CWarnerValveControllerDeviceNet, 659
IsDeviceHighSpeed	IsValveOpenInAnalogModeEvent
CMcsUsbNet, 329	CWarnerValveControllerDeviceNet, 671
IsDeviceHighSpeedCapable	IsValveOpenInDigitalMode
CMcsUsbNet, 329	CWarnerValveControllerDeviceNet, 659
IsDeviceTypeOf	IsValveOpenInDigitalModeEvent
CMcsUsbListNet, 317	CWarnerValveControllerDeviceNet, 671
IsDigitalChannelDedicated	
CMcsUsbDacqNet::CHWInfo, 188	Kelvin
IsDigitalOutPortInverted	Mcs::Usb, 53
CWarnerValveControllerDeviceNet, 658	
IsDigitalOutPortInvertedEvent	LastPosition
CWarnerValveControllerDeviceNet, 670	Mcs::Usb, 57, 73–75, 84, 86, 93
IsEqual	LegacyHasOutputMap
CFilterCoefficientsNet, 136	CStg200xBasicNet, 542
IsExceptionsEnabled	LegacyMeaUsb
CMcsUsbNet, 330	Mcs::Usb, 77
IsGateFloating	LIH30_ADC_Channel_EnumNet
CCMOSMea_FunctionNet, 111	Mcs::Usb, 68
IsHeadstageAvailable	LIH30_DAC_Channel_EnumNet
CSCUFunctionNet, 514	Mcs::Usb, 69
IsHeadstageAvailableEvent	LIH30_EPC10_Bus_EnumNet
CSCUFunctionNet, 519	Mcs::Usb, 69
IsHighCurrentMode	LIH30ADCCtrl
CWarnerUssingFunctionNet, 642	Mcs::Usb, 67, 68
IsHSPowered	LIH30ADCModulesGroup
CSCUFunctionNet, 515	Mcs::Usb, 61
IsInDacqLegacyMode	LIH30Interfaceboard
CSCUFunctionNet, 515	Mcs::Usb, 67
IsInternalCalibrationFinished	LIH30TestADCGroup
CTEERFunctionNet, 607	Mcs::Usb, 61
	LIH30UserADCGroup
CWarnerUssingFunctionNet, 642	Mcs::Usb, 61
IsPlateTypeValid	ListModeSendStart
CMultiwellDeviceNet, 388	CStg200xBasicNet, 542
IsPulseEnabled	ListModeSendStop
CWarnerUssingFunctionNet, 642	CStg200xBasicNet, 543
IsPumpMotorOn	ListOfChangedTriggers
CRoboFluidDeviceNet, 489	StgStatusNet, 712
IsQueueEnabled	LoadPressure
CRoboDeviceNet, 478	CPPCFunctionNet, 431
IsQueueStarted	LoadUserFirmware
CRoboDeviceNet, 478	CMcsUsbFactoryNet, 304, 305
IsRunning	LoadValveTable
CMeaCleanDeviceNet, 346	CWarnerValveControllerDeviceNet, 659
CMeaCoatDeviceNet, 351	Lock
IsSamplingFinished	Mcs::Usb, 76
CTEERFunctionNet, 607	LockPlateClamp
IsUserTriggerEnabled	CMultiwellDeviceNet, 389
CLIH3DeviceNet, 200	
IsValveDigitalInInverted	Lowpass Mos: Usb. 66
CWarnerValveControllerDeviceNet, 658	Mcs::Usb, 66
IsValveDigitalInInvertedEvent	LowSpeed
CWarnerValveControllerDeviceNet, 670	Mcs::Usb, 70
IsValveOpen	M4M Dock
CWarnerValveControllerDeviceNet, 659	-
	Mcs::Usb, 67

M4M_ReUse	ALTERA_TARGET3, 56
Mcs::Usb, 67	AlwaysOn, 79
m_Bottom	Ampere, 53, 91
CCMOSMeaDeviceNet::CRegionOfInterestRect,	AMS_Dongle, 77
446	Analog, 76, 94
m_Left	AnalogGroup, 61
CCMOSMeaDeviceNet::CRegionOfInterestRect,	AnalogOut_DAC_Range_EnumNet, 52
446	AnalogSource_HS1, 53
m_pGilsonDevice	AnalogSource_HS2, 53
CGilsonDeviceNet, 169	AnalogSource IF, 53
,	- '
m_pMcsBus_MotorControlNet	AnalogSourceEnumNet, 53
CRoboFluidDeviceNet, 490	AnalogUnitEnumNet, 53
m_pMcsUsb	Any, 77, 91
CMcsUsbFunctionNet, 309	Armed, 85
m_pMcsUsbFunction	ASMedia, 90
CMcsUsbFunctionNet, 309	AudioTestChannelGroup, 61, 71, 82, 93
m_pRoboFluidDevice	Aux, 63
CRoboFluidDeviceNet, 490	AuxIn, 57, 71, 73, 74, 83, 86, 93
m_Right	AuxPort, 59
CCMOSMeaDeviceNet::CRegionOfInterestRect,	B, 69
446	Bessel, 66
m Top	BMI, 87
CCMOSMeaDeviceNet::CRegionOfInterestRect,	Bootstrap, 55
446	BootstrapOtherCypress, 55
Manual	Both, 91
Mcs::Usb, 76, 84, 94	Break, 81
Manufacturer	BUSOMCSBUSO, 54
CMcsUsbListEntryNet, 314	BUS0MCSBUS1, 54
MaxBitNumber	BUS0MCSBUS10, 54
DigitalSource< digitalsourceenum >, 682	BUS0MCSBUS11, 54
DigitalSourceGeneral, 683, 684	BUS0MCSBUS12, 54
MaxBitNumberStatic	BUS0MCSBUS13, 54
DigitalSource< digitalsourceenum >, 682	BUS0MCSBUS14, 54
MBC08	BUS0MCSBUS15, 54
Mcs::Usb, 78	BUS0MCSBUS2, 54
MbcChannelStateEnumNet	BUS0MCSBUS3, 54
Mcs::Usb, 69	BUS0MCSBUS4, 54
MbcChargingModeEnumNet	BUS0MCSBUS5, 54
Mcs::Usb, 70	BUS0MCSBUS6, 54
MbcRatedCapacityEnumNet	BUS0MCSBUS7, 54
Mcs::Usb, 70	
	BUSOMCSBUS8, 54
MC_Card	BUSOMCSBUS9, 54
Mcs::Usb, 77	BUS1MCSBUS0, 54
MCS	BUS1MCSBUS1, 54
Mcs::Usb, 91	BUS1MCSBUS10, 54
Mcs, 22	BUS1MCSBUS11, 54
Mcs::Usb, 23	BUS1MCSBUS12, 54
A, 69	BUS1MCSBUS13, 54
AccelOnly, 91	BUS1MCSBUS14, 54
AdapterTypeEnumNet, 52	BUS1MCSBUS15, 54
ALA_VC3, 77, 91	BUS1MCSBUS2, 54
ALA_VC3_DEVICE, 63	BUS1MCSBUS3, 54
ALTERA, 55	BUS1MCSBUS4, 54
ALTERA_BASE, 56	BUS1MCSBUS5, 54
ALTERA_BOOTSTRAP, 56	BUS1MCSBUS6, 54
ALTERA_GOLD, 55	BUS1MCSBUS7, 54
ALTERA_TARGET1, 56	BUS1MCSBUS8, 54
ALTERA_TARGET2, 56	BUS1MCSBUS9, 54

BUS2MCSBUS0, 55	CMOSMeaPacketFrameContextGroupEnumNet,
BUS2MCSBUS1, 54	59
BUS2MCSBUS10, 55	CMOSMeaSTG1DACSignalEnumNet, 60
BUS2MCSBUS11, 55	CMOSMeaValueUnitEnumNet, 60
BUS2MCSBUS12, 55	CommaPositionA, 65
BUS2MCSBUS13, 55	CommaPositionB, 65
BUS2MCSBUS14, 55	csCapacityTestDischarge, 69
BUS2MCSBUS15, 55	csCapacityTestPrecharge, 69
BUS2MCSBUS2, 54	csCharge, 69
BUS2MCSBUS3, 54	csDischarge, 69
BUS2MCSBUS4, 54	csError, 69
BUS2MCSBUS5, 55	csldleChargeFinished, 69
BUS2MCSBUS6, 55	csldleNoBattery, 69
BUS2MCSBUS7, 55	csRefreshBattery, 69
BUS2MCSBUS8, 55	CurrentClamp, 90
BUS2MCSBUS9, 55	CurrentMeasure, 57
BUSNUMBER0, 54	Cypress, 91
BUSNUMBER1, 54	Cypress FX1, 77
BUSNUMBER2, 54	Cypress_FX2, 77
Butterworth, 66	Cypress FX3, 77
Campden_Ci4600EphysVideoDataIntegrator, 77	DAC1Channel, 60
CatchAmp, 75	DAC2Channel, 60
CFirmwareDestinationNet, 53	DAC3Channel, 60
channeldata_current, 85	DAC4Channel, 60
channeldata_current_always_boost, 86	DACQ1DigitalGroup, 61
channeldata_current_always_boost_own_sync, 86	DacqGroupChannelEnumNet, 60
channeldata_current_own_boost_gnd_sync, 86	DacqMeaGroupTypeEnumNet, 61
channeldata_current_own_sync, 86	DacqTrigger, 64
channeldata_positive_current, 86	DataModeEnumNet, 61
channeldata_positive_current_own_boost_gnd_sync,	DEST_FX3_TARGET_MASK, 56
86	DEST_TARGET1, 56
channeldata_positive_current_own_sync, 86	DEST_TARGET10, 56
channeldata positive voltage, 86	DEST_TARGET11, 56
channeldata_voltage, 85	DEST_TARGET12, 56
ChannelPIC, 55	DEST_TARGET13, 56
ChannelTest, 78	DEST_TARGET14, 56
ChecksumAndPacketCounter, 60	DEST_TARGET15, 56
Ci4600Intan, 52	DEST_TARGET2, 56
ClampModeCurrent, 87	DEST_TARGET3, 56
ClampModeInternalCalibration, 87	DEST_TARGET4, 56
ClampModeOpen, 87	DEST_TARGET5, 56
ClampModeVoltage, 87	DEST_TARGET6, 56
Close, 76	DEST_TARGET7, 56
CmosMea, 68	DEST_TARGET8, 56
CMOSMEA5000DigitalSourceEnumNet, 57	DEST TARGET9, 56
CMOSMeaBathModeEnumNet, 57	DEST_TARGET_MASK, 56
CmosMeaHeadstage, 67	DEVICE_NOT_FOUND, 65
CMOSMeaHeadstage1NCBathCurrentEnumNet,	DeviceEnumNet, 61
57	DeviceHasNoHeadstage, 67, 68, 82
CMOSMeaHeadstage1NCCol2CurrentEnumNet,	DeviceNotConnected, 66, 67, 81
58	DeviceRunStatus, 57, 71, 74, 75, 83, 86, 93
CMOSMeaHeadstage1NChipTempEnumNet, 58	DigDataFromReceiver, 93
CMOSMeaHS1SidebandEnumNet, 58	Digital, 76, 94
CMOSMeaHS1TriggerStatusEnumNet, 58	DigitalData, 57, 71, 74, 75, 83, 86, 93
CmosmealFB2, 67	DigitalDatastreamEnableEnumNet, 63
CMOSMealFDigChannelEnumNet, 59	DigitalGroup, 61
CMOSMeaInterfaceADCEnumNet, 59	DigitalIn, 57, 63, 71, 73, 74, 83, 86, 93
CmosMeaInterfaceboard, 67	DigitalInOfOutPort, 57, 71, 73, 74, 83, 86, 93
	J : - : - : - : - : - : - : - :

DigitalInPort, 59	FPGA10_GOLD, 56
DigitalInReserverd, 63	FPGA11, 55
DigitalMux, 59	FPGA11_BASE, 56
DigitalOut, 63	FPGA11_GOLD, 56
DigitalOutReg, 59	FPGA12, 55
DigitalOutReserved, 63	FPGA12_BASE, 56
DigitalOutStimulator, 57, 71, 74, 75, 83, 86, 93	FPGA12 GOLD, 56
DigitalPulse, 57, 71, 73, 74, 83, 86, 93	FPGA13, 55
DigitalReg, 59	FPGA13_BASE, 56
<u> </u>	
DigitalStimulatorTriggerEventEnumNet, 63	FPGA13_GOLD, 56
DigitalStimulatorTriggerSlopeEnumNet, 64	FPGA14, 55
DigitalTargetEnumNet, 64	FPGA14_BASE, 56
Digout, 64	FPGA14_GOLD, 56
DigOutStim, 63	FPGA15, 55
DigOutStimulatorStartTrigger, 64	FPGA15_BASE, 56
DigOutStimulatorStopTrigger, 64	FPGA15_GOLD, 56
Digstream, 64	FPGA16, 55
DigStreamFromReceiver, 93	FPGA16_BASE, 56
DigStreamToReceiver, 64	FPGA16_GOLD, 56
Dilutor, 79	FPGA2, 55
DongleS, 78	FPGA2_BASE, 56
Dotriapot, 77	FPGA2_GOLD, 55
DoubleToInt, 66	FPGA3, 55
DownloadOnly, 85	FPGA3_BASE, 56
DSP, 53	FPGA3_GOLD, 55
DSPAnalogGroup, 71, 82	FPGA4, 55
DSPDataGroup, 60, 92	FPGA4_BASE, 56
DSPDigitalGroup, 71, 82	FPGA4 GOLD, 55
eCube, 78, 82	FPGA5, 55
eCubeHeadstage, 67	FPGA5 BASE, 56
ElectrodeDacMuxEnumNet, 64	FPGA5_GOLD, 55
ElectrodeModeEnumNet, 64	FPGA6, 55
ElectrodeOffset, 90	FPGA6_BASE, 56
emAutomatic, 65	FPGA6_GOLD, 55
emManual, 65	FPGA7, 55
Encapsulator, 79	FPGA7_BASE, 56
enCMosMeaChipType, 65	FPGA7_GOLD, 56
EnSTG200x_STATUS, 65	FPGA8, 55
EOFAndCRC, 60	FPGA8_BASE, 56
ExtendedMultiFileMode, 84	FPGA8_GOLD, 56
ExternBCTester, 78	FPGA9, 55
ExternDTester, 78	FPGA9_BASE, 56
ExternSTester, 78	FPGA9_GOLD, 56
Falling, 64	FPGA_BASE, 56
FCB, 77	FPGA BOOTSTRAP, 56
FCX, 77	FPGA_GOLD, 55
Feedback, 57, 71, 73, 74, 83, 86, 93	FPGA_NORMAL, 53
FeedbackHigh, 63	FpgaldEnumNet, 66
FeedbackLow, 63	FrameContextGroup, 61
FeedbackReg, 59	FullCharge, 70
FilterAttributeEnumNet, 65	FullSpeed, 70
FilterBandEnumNet, 65	FunkDongleS, 78
FilterCalculationDirectionEnumNet, 66	_
	Gain100, 85
FilterFamilyEnumNet, 66	Gain1000, 85
FilterTypeEnumNet, 66	Gain200, 85
Finished, 85	Gain500, 85
FPGA10, 55	Gate, 80
FPGA10_BASE, 56	Gated_High_Active, 79

Gated_Low_Active, 79	Hs1Trigger, 63
GE2100, 78	HS1Trigger10Status, 72
Generic, 77	HS1Trigger11Status, 72
Graphene_ASIC, 78	HS1Trigger12Status, 72
GrapheneASIC, 82	HS1Trigger13Status, 72
GrapheneASICHeadstage, 67	HS1Trigger14Status, 72
GrapheneProjectTestDevice, 79	HS1Trigger15Status, 72
Ground, 57, 64	HS1Trigger16Status, 72
GyroOnly, 91	HS1Trigger17Status, 72
Hardware, 66	HS1Trigger18Status, 72
Headstage1NCBathCurrentGroup, 60	HS1Trigger1Status, 57, 71, 73, 74
Headstage1NCCol2CurrentGroup, 60	HS1Trigger2Status, 57, 71, 73, 74
Headstage1NChipTempGroup, 61	HS1Trigger3Status, 57, 71, 73
HeadstageElectrodeGroup, 60	HS1Trigger4Status, 71, 73
HeadstageldEnumNet, 67	HS1Trigger5Status, 71, 73
HEKA_LIH3_DEVICE, 62	HS1Trigger6Status, 71, 73
HekaEPC10Double, 77	HS1Trigger7Status, 71
HekaEPC10Quadro, 77	HS1Trigger8Status, 72
HekaEPC10Single, 77	HS1Trigger9Status, 72
HekaEPC10Triple, 77	Hs2Digital, 63
HekaEPCLite, 77	HS2DigitalData1, 86
HekalTEV100, 77	HS2ElectrodeGroup, 71
HekaLIH30, 77	HS2Sideband1, 72, 74, 75
HekaLIH406, 77	HS2Sideband10, 73
HekaLIH816, 77	HS2Sideband11, 73
HekaPG610, 77	HS2Sideband12, 73
HekaPG611, 77	HS2Sideband13, 73
HekaPG612, 77	HS2Sideband14, 73
HekaPG618, 77	HS2Sideband15, 73
HekaPG690, 77	HS2Sideband16, 73
HiClamp, 79	HS2Sideband17, 73
HiClamp4Uart, 79	
•	HS2Sideband18, 73 HS2Sideband2, 73–75
Highpass, 66 HighSpeed, 70	
•	HS2Sideband3, 73, 74
HLA, 77 Hs1Digital, 63	HS2Sideband4, 73, 74 HS2Sideband5, 73, 74
HS1DigitalData1, 86	HS2Sideband6, 73, 74
HS1ElectrodeGroup, 71	• •
HS1Sideband1, 57, 72–74	HS2Sideband7, 73
	HS2Sideband8, 73
HS1Sideband10, 72 HS1Sideband11, 72	HS2Sideband9, 73 Hs2SidebandHigh, 63
HS1Sideband12, 72	Hs2SidebandLow, 63
HS1Sideband13, 72	Hs2Trigger, 63
HS1Sideband14, 72	HS2Trigger10Status, 72
HS1Sideband15, 72	HS2Trigger11Status, 72
HS1Sideband16, 72 HS1Sideband17, 72	HS2Trigger12Status, 72
•	HS2Trigger13Status, 72
HS1Sideband18, 72	HS2Trigger14Status, 72
HS1Sideband2, 57, 72–74	HS2Trigger15Status, 72
HS1Sideband3, 57, 72, 73	HS2Trigger16Status, 72
HS1Sideband4, 72, 73	HS2Trigger12Status, 72
HS1Sideband5, 72, 73	HS2Trigger18Status, 72
HS1Sideband6, 72, 73	HS2Trigger1Status, 72–74
HS1Sideband7, 72	HS2Trigger2Status, 72–74
HS1Sideband8, 72	HS2Trigger3Status, 72, 73
HS1Sideband9, 72	HS2Trigger4Status, 72, 74
Hs1SidebandHigh, 63	HS2Trigger5Status, 72, 74
Hs1SidebandLow, 63	HS2Trigger6Status, 72, 74

HS2Trigger7Status, 72	MCS_ENCAPSULATOR_DEVICE, 62
HS2Trigger8Status, 72	MCS_EXTERN_BC_TESTER_DEVICE, 62
HS2Trigger9Status, 72	MCS_EXTERN_D_TESTER_DEVICE, 62
Idle, 85	MCS_FCX_DEVICE, 62
IFB2GoldenInterfaceboard, 67	MCS_FYI_DEVICE, 62
IFB30GoldenInterfaceboard, 67	MCS_GENERIC_DEVELOPMENT_DEVICE, 62
IFChannel1, 59	MCS_HICLAMP_DEVICE, 62
IFChannel2, 59	MCS_HLA_DEVICE, 62
IFChannel3, 59	MCS_MBC08_DEVICE, 62
IFChannel4, 59	MCS_MC_STIMULUS_DEVICE, 62
IFChannel5, 59	MCS_MCCARD_DEVICE, 62
IFChannel6, 59	MCS_MEA_CLEAN_DEVICE, 62
IFChannel7, 59	MCS_MEA_COAT_DEVICE, 62
IFChannel8, 59	MCS_MEA_DEVICE, 62
IFDigChannelsGroup, 61, 71, 82, 93	MCS_MEA_IMPEDANCE_DEVICE, 62
Ignore, 80	MCS_MEA_SWITCH_DEVICE, 62
IM16KRC, 79	MCS_MEASURETABLE_DEVICE, 62
IM16S16KRA, 79	MCS_MEAUSB_DEVICE, 62
IM16S8KRA, 79	MCS_NF_GEN_DEVICE, 62
IM64KRB, 79	MCS_OCTOPOT_DEVICE, 62
IM64KRC, 79	MCS_OKUVISION_STIMULATOR_DEVICE, 62
Input, 76	MCS_PATCHSERVER_DEVICE, 62
Intel, 90	MCS_PATHIDENT_DEVICE, 62
InterfaceADCGroup, 60, 71, 82, 91 InterfaceBoard2, 67	MCS_PCI_BUS, 70
IntToDouble, 66	MCS_PCX_DEVICE, 62 MCS_PEDOTER_DEVICE, 62
InvitroSignalCollectorUnit, 68	MCS_PERISTALTIC_PUMP_DEVICE, 62
InvivoSignalCollectorUnit, 68	MCS_PGA_DEVICE, 62
IoVoltageEnumNet, 68	MCS_PPC_DEVICE, 62
IS32KRA, 79	MCS_PPS5_DEVICE, 62
Kelvin, 53	MCS_PPS_DEVICE, 62
LastPosition, 57, 73–75, 84, 86, 93	MCS_RETINA_AMS_DONGLE, 62
LegacyMeaUsb, 77	MCS_RETINA_LED_DEVICE, 62
LIH30 ADC Channel EnumNet, 68	MCS_ROBO_DEVICE, 62
LIH30_DAC_Channel_EnumNet, 69	MCS_ROBOINJECT_DEVICE, 62
LIH30_EPC10_Bus_EnumNet, 69	MCS_ROBOOCYTE2_DEVICE, 62
LIH30ADCCtrl, 67, 68	MCS_SAFEIS_DEVICE, 62
LIH30ADCModulesGroup, 61	MCS_SMARTIMPLANT_DEVICE, 62
LIH30Interfaceboard, 67	MCS SOFTWARE DONGLE DEVICE, 62
LIH30TestADCGroup, 61	MCS STG DEVICE, 62
LIH30UserADCGroup, 61	MCS_SW2TO64_DEVICE, 62
Lock, 76	MCS_TCX_DEVICE, 62
Lowpass, 66	MCS_TERSENS_DEVICE, 62
LowSpeed, 70	MCS_UNDEFINED_BUS, 70
M4M_Dock, 67	MCS_USB_BUS, 70
M4M_ReUse, 67	MCSBUS0, 54
Manual, 76, 84, 94	MCSBUS1, 53
MBC08, 78	MCSBUS10, 54
MbcChannelStateEnumNet, 69	MCSBUS11, 54
MbcChargingModeEnumNet, 70	MCSBUS12, 54
MbcRatedCapacityEnumNet, 70	MCSBUS13, 54
MC_Card, 77	MCSBUS14, 54
MCS, 91	MCSBUS15, 54
MCS_ANY_BUS, 70	MCSBUS2, 53
MCS_CHANNELTEST_DEVICE, 62	MCSBUS3, 53
MCS_DEVICE_ANY, 62	MCSBUS4, 53
MCS_DEVICE_USB, 62	MCSBUS5, 53
MCS_DEVICE_USB_CYPRESS, 63	MCSBUS6, 54

M000U07 54	MEAGO 50
MCSBUS7, 54	MEA60, 52
MCSBUS8, 54	MEA_2_252_2, 52
MCSBUS9, 54	MEA_2_252_2_6Well, 52
McsBusTypeEnumNet, 70	MEA_2_252_2_9Well, 52
McsUsbSpeedEnumNet, 70	MEA_2_252_2_Test, 52
MCU1, 53	MEA_Clean, 78
ME128, 78	MEA_Coat, 78
ME16, 78	MEA_Impedance, 78
ME2100, 78	MEA Sanofi, 78
Me2100 32PICiCE40, 81	MEA Switch, 78
Me2100 32PICiCE40Headstage, 67	MEA Switch 2 1, 78
Me2100_32Xilinx, 81	MEA_Switch_4_2, 78
Me2100_32XilinxHeadstage, 67	MeaLayoutEnumNet, 74
Me2100Graphene16_32, 82	MeasureTable, 79
Me2100Graphene16_32Headstage, 67	MicroAmpere, 60
. –	•
Me2100Interfaceboard, 67	MilliDegreeCelsius, 60
Me2100InvitroSignalCollectorUnit, 67	mIMEA60, 74
Me2100InvivoSignalCollectorUnit, 67	mlUnknown, 74
Me2100UPA32, 81	ModulA_ADC0, 68
Me2100UPA32Headstage, 67	ModulA_ADC1, 68
ME256, 78	ModulA_ADC2, 68
ME32, 78	ModulA_ADC3, 68
ME64, 78	ModulA_DAC0, 69
MEA1060, 78	ModulA_DAC1, 69
MEA120, 52	ModulB_ADC0, 68
MEA2100, 78	ModulB_ADC1, 68
Mea2100, 67	ModulB_ADC2, 68
MEA2100_256, 78	ModulB ADC3, 68
Mea2100_256, 68	ModulB DAC0, 69
MEA2100_256DacqGroupChannelEnumNet, 71	ModulB DAC1, 69
MEA2100_256DigitalSourceEnumNet, 71	ModulC_ADC0, 68
Mea2100_256Headstage, 67	ModulC_ADC1, 68
Mea2100_256Interfaceboard, 67	ModulC_ADC2, 68
	ModulC_ADC3, 68
MEA2100_32, 78	
MEA2100_Lite, 78	ModulC_DAC0, 69
Mea2100_Lite, 68	ModulC_DAC1, 69
MEA2100_Mini, 78	ModulD_ADC0, 69
MEA2100_Mini_Usb_develop, 78	ModuID_ADC1, 69
MEA2100BetaScreen, 78	ModulD_ADC2, 69
Mea2100BetaScreen, 81	ModuID_ADC3, 69
Mea2100BetaScreenHeadstage, 67	ModuID_DAC0, 69
MEA2100DigitalSourceEnumNet, 73	ModulD_DAC1, 69
Mea2100Headstage, 66	Monitor, 84
Mea2100Interfaceboard, 66	Movement, 81
Mea2100LiteHeadstage, 67	Multiboot, 79
Mea2100Mini120, 82	MultiFileMode, 84
Mea2100Mini120Headstage, 67	Multiwell, 68, 78
Mea2100Mini60ECP5, 82	Multiwell96, 52
Mea2100Mini60ECP5Headstage, 67	Multiwell ICC, 78
Mea2100Mini60PICiCE40, 81	Multiwell_MEA_Mini, 78
Mea2100Mini60PICiCE40Headstage, 67	MultiwellDigitalSourceEnumNet, 74
Mea2100MultiwellIFB2, 67	MultiwellHeadstage, 66
Mea2100STG, 66	MultiwellInterfaceboard, 66
MEA252, 52	MultiwellMini, 82
MEA256, 78	MultiwellOntoStim, 77
MEA2x32, 52	MultiwellOptoStim, 77
MEA2x60, 52	MultiwellPlateTypeEnumNet, 75
MEA32, 52	Mux, 63

	D100
MuxOtherDevice, 63	PIC3, 55
Nanion, 87	PIC4, 55
NanoAmpere, 60	PIC5, 55
NanoVolt, 60	PIC6, 55
NCBathCurrent, 58	PIC7, 55
NCCol2Current, 58	PIC8, 55
NChipTemperature, 58	PIC9, 55
Neptun, 79	PicoAmpere, 60
NeuroChip, 79	Plate_24W030MGA, 75
NeurochipConfig, 79	Plate_24W300_30_1152GBA, 75
NF_Gen, 79	Plate_24W300_30GBA, 75
NineWell, 86	Plate_24W300_30GBB, 75
nMos16LV, 65	Plate_24W300_30GMA, 75
nMos32LV, 65	Plate_24W700_100FMA, 75
nMos36LN, 65	Plate_24W700_100FMB, 75
nMos64LN, 65	Plate_24W700_100FMC, 75
No_Plate, 75 None, 52, 63, 77, 85, 90, 91	Plate_24W700_100PBA, 75 Plate_72W500_100FMA, 75
Normal, 75	Plate 72W500_100PMA, 75
NOT CONNECTED, 65	Plate 96W300 80 1152FMA, 75
NotApplicable, 52	Plate 96W400 80 1152FMB, 75
NoUnit, 60	Plate_96W700_100FMA, 75
NTC10K, 87	Plate_96W700_100FMB, 75
Octopot, 77	Plate_96W700_100GBA, 75
Off, 79, 81, 91	Plate_96W700_100GBB, 75
off, 80	Plate_96W700_100GBC, 75
OK, 65	Plate_96W700_100GBD, 75
Okuvision_Stimulator, 78	Plate_96W700_100GMA, 75
OnChannelData, 94	Plate Dummy, 75
OnDeviceArrivalRemoval, 94	Plate_Dummy_126, 75
One, 57, 71, 73, 74, 83, 86, 93	PlateClampEnumNet, 76
OnError, 94	PlateClampLockEnumNet, 76
OnMcsUsbDeviceState, 94	PlusMinus10Volts, 53
OnMcsUsbDeviceStateCallback, 95	PlusMinus2Comma5Volts, 53
OnMwPollStatus, 95	PlusMinus5Volts, 53
OnStg200xDataHandler, 95	PortDirectionEnumNet, 76
OnStg200xErrorHandler, 95	Pos900, 79
OnStgPollStatus, 95	PositionBase, 79
OnUpdateFirmwareProgress, 95	PositionIIBase, 79
OnUpdateFirmwareStatusChange, 95	PositionIICentralUnit, 79
Open, 76	PositionImp, 79
OpenClamp, 90	PostCommaA, 65
Output, 76	PostCommaB, 65
PacketFrameContextGroup, 61, 71, 82, 93	PP_Pump_Mode_Type_EnumNet, 76
PatchServAdcModeEnumNet, 75	PPC, 78
PatchServer, 79	PPS2, 78
PathIdent, 79	PPS4plus1, 78
PC, 94	PPS5, 78
PCI, 91	PPS5_DIG, 78
PCX, 77	PreCommaA, 65
PeriodicPulse, 63	PreCommaB, 65
PeristalticPump, 78	ProductIdEnumNet, 76
PGA, 77	PT100, 87
PIC, 55	PT1000, 87
PIC10, 55	PulseGenerator, 57, 71, 74, 75, 83, 86, 93
PIC11, 55	PulseGenerator_Mode_EnumNet, 79
PIC12, 55	rawdata, 86
PIC2, 55	RC, 66

rc100mAh, 70	SCU1HS2Sideband2, 83
rc200mAh, 70	SCU1HS2Sideband2TriggerStatus, 83
rc300mAh, 70	SCU1HS3Sideband1, 83
rc30mAh, 70	SCU1HS3Sideband1TriggerStatus, 83
rcGreater300mAh, 70	SCU1HS3Sideband2, 83
Rectangle, 88	SCU1HS3Sideband2TriggerStatus, 83
Ref16, 80	SCU1HS4Sideband1, 83
Ref24, 80	SCU1HS4Sideband1TriggerStatus, 83
Ref32, 80	SCU1HS4Sideband2, 83
Ref8, 80	SCU1HS4Sideband2TriggerStatus, 83
Reference, 81	SCU1LEDStim1, 83
ReferenceElectrodeModeEnumNet, 80	SCU1LEDStim1TriggerStatus, 83
ReferenceElectrodeSwitchPositionEnumNet, 80	SCU1LEDStim2, 83
RegisterHigh, 63	SCU1LEDStim2TriggerStatus, 83
RegisterLow, 63	SCU1LEDStim3, 83
Regular, 87	SCU1LEDStim3TriggerStatus, 83
Renesas, 90	SCU1LEDStim4, 83
Reserved1, 87	SCU1LEDStim4TriggerStatus, 83
Reserved2, 87	SCU2ElectrodeGroupHS1, 82
Reserved3, 87	SCU2ElectrodeGroupHS2, 82
Reserved4, 87	SCU2ElectrodeGroupHS3, 82
Reserved5, 87	SCU2ElectrodeGroupHS4, 82
Restart, 80	SCU2HS1Sideband1, 83
Retina LED, 77	SCU2HS1Sideband1TriggerStatus, 83
RetriggerActionEnumNet, 80	
Rising, 64	SCU2HS1Sideband2TriggerStatus 82
•	SCU2HS1Sideband2TriggerStatus, 83
RoboCurrentModeEnumNet, 80	SCU2HS2Sideband1, 84
Robolnject, 79	SCU2HS2Sideband1TriggerStatus, 83
Robooycte2, 79	SCU2HS2Sideband2, 84
RoboStatusEventDelegate, 95	SCU2HS2Sideband2TriggerStatus, 83
Running, 85	SCU2HS3Sideband1, 84
SafelS, 79	SCU2HS3Sideband1TriggerStatus, 83
SampleDstSize16, 81	SCU2HS3Sideband2, 84
SampleDstSize32, 81	SCU2HS3Sideband2TriggerStatus, 83
SampleDstSizeNet, 81	SCU2HS4Sideband1, 84
SampleSize16Signed, 81	SCU2HS4Sideband1TriggerStatus, 83
SampleSize16Unsigned, 81	SCU2HS4Sideband2, 84
SampleSize24Signed, 81	SCU2HS4Sideband2TriggerStatus, 83
SampleSize24Unsigned, 81	SCU2LEDStim1, 84
SampleSize32Signed, 81	SCU2LEDStim1TriggerStatus, 83
SampleSize32Unsigned, 81	SCU2LEDStim2, 84
SampleSize64Signed, 81	SCU2LEDStim2TriggerStatus, 83
SampleSize64Unsigned, 81	SCU2LEDStim3, 84
SampleSizeNet, 81	SCU2LEDStim3TriggerStatus, 83
SBSVector1, 58	SCU2LEDStim4, 84
SBSVector2, 58	SCU2LEDStim4TriggerStatus, 83
SBSVector3, 58	SCU_HeadstageIdEnumNet, 81
SBSVector4, 58	SCUDacqGroupChannelEnumNet, 82
SCU1ElectrodeGroupHS1, 82	SCUDigitalSourceEnumNet, 82
SCU1ElectrodeGroupHS2, 82	Settings, 94
SCU1ElectrodeGroupHS3, 82	Signed_16bit, 61
SCU1ElectrodeGroupHS4, 82	Signed_24bit, 61
SCU1HS1Sideband1, 83	Signed_32bit, 61
SCU1HS1Sideband1TriggerStatus, 83	Sine, 88
SCU1HS1Sideband2, 83	Single, 80
SCU1HS1Sideband2TriggerStatus, 83	SingleWell, 86
SCU1HS2Sideband1, 83	SixWell, 86
SCU1HS2Sideband1TriggerStatus, 83	SmartImplant, 79

SOFAndCTRLword, 60	TBSI_5, 52
Software, 66	TBSI_63, 52
SoftwareDongle, 79	TBSI_Dacq, 78
Standby, 81, 90	TBSI_DACQDigitalSourceEnumNet, 86
Start, 63	TBSI_Reserved, 52
State, 91	TbsiDacq, 68
STG, 77	TbsiDacqHeadstage, 67
Stg1, 64	TbsiDacqInterfaceboard, 67
STG1DACSignalGroup, 61, 71, 82	TC01, 77
STG1SidebandsGroup, 61, 71, 82	TC02, 77
STG1TriggerStatusGroup, 61, 71, 82	TCX, 77
Stg2, 64	TcxDeviceTypeEnumNet, 87
Stg200xDigoutModeEnumNet, 84	TcxSensorTypeEnumNet, 87
Stg200xMultiFileSubmodeEnumNet, 84	TeerClampModeEnumNet, 87
Stg200xSegmentFlagsEnumNet, 84	TeerWaveformEnumNet, 87
Stg200xTriggerStatusEnumNet, 85	Tersens, 77
STG2DACSignalGroup, 71, 82	Test_ADC_EPC10, 68
STG2SidebandsGroup, 71, 82	Test_DAC_EPC10, 69
STG2TriggerStatusGroup, 71, 82	Timestamp, 60
Stg3, 64	TouchTest, 94
STG3008_FA, 77	Triggerbox_AMS, 78
Stg3008FilterAmpAmplificationEnumNet, 85	Triggerbox AMS3, 78
STG4002, 77	Triggerbox IMS, 78
STG4002_opto, 77	Triggerbox_R5, 78
STG4004, 77	TriggerOnly, 85
STG4004_opto, 77	TriggerSourceEnumNet, 88
STG4008, 77	TriggerStatus1, 59
STG4008_opto, 77	TriggerStatus2, 59
STG400x, 77	TriggerStatus3, 59
STG400x_opto, 77	TriggerStatus4, 59
STG5, 77	tsAuxIn1, 89
STG_DestinationEnumNet, 85	tsAuxIn2, 89
StgListModeTrigger, 64	tsDACQCy1Dev1Runs, 90
StgTrigger, 64	tsDACQCy1Dev2Runs, 90
Stimulation, 57	tsDACQCy2Dev1Runs, 90
StimulationLayoutConfigurationEnumNet, 86	tsDACQCy2Dev2Runs, 90
Stop, 63, 76, 80	tsDigitalIn1, 88
StorageCharge, 70	tsDigitalIn10, 88
SubtractFromAll, 80	tsDigitalIn11, 88
SubtractFromAllOther, 80	tsDigitalIn12, 88
SubtractFromReferenceElectrodeOnly, 80	tsDigitalIn13, 88
SubtractionOff, 80	tsDigitalIn14, 88
SuperSpeed, 70	tsDigitalIn15, 88
Sw2to64, 78	tsDigitalIn16, 88
SYNCOUT1, 84	tsDigitalIn17, 88
SYNCOUT2, 84	tsDigitalIn18, 88
SYNCOUT3, 84	tsDigitalIn19, 88
SYNCOUT4, 84	tsDigitaIIn2, 88
SYNCOUT5, 84	tsDigitalIn20, 88
SYNCOUT6, 84	tsDigitalIn21, 88
SYNCOUT7, 84	tsDigitalIn22, 88
SYNCOUT8, 84	tsDigitalIn23, 88
syncoutdata, 85	tsDigitalIn24, 88
SyncStart, 85	tsDigitalIn25, 88
Table, 94	tsDigitalIn26, 88
TBSI_127, 52	tsDigitalIn27, 88
TBSI_15, 52	tsDigitalIn28, 88
TBSI_31, 52	tsDigitaIIn29, 88

tsDigitalIn3, 88	tsFeedback24, 89
tsDigitalIn30, 88	tsFeedback25, 89
tsDigitalIn31, 88	tsFeedback26, 89
tsDigitalIn32, 88	tsFeedback27, 89
tsDigitalIn4, 88	tsFeedback28, 89
tsDigitalIn5, 88	tsFeedback29, 89
tsDigitalIn6, 88	tsFeedback3, 88
tsDigitalIn7, 88	tsFeedback30, 89
tsDigitalIn8, 88	tsFeedback31, 89
tsDigitalIn9, 88	tsFeedback32, 89
tsDigitalPuse0, 89	tsFeedback4, 89
tsDigitalPuse1, 89	tsFeedback5, 89
tsDigitalPuse10, 89	tsFeedback6, 89
tsDigitalPuse11, 89	tsFeedback7, 89
tsDigitalPuse12, 89	tsFeedback8, 89
tsDigitalPuse13, 89	tsFeedback9, 89
_	tsNone, 88
tsDigitalPuse15, 89	
tsDigitalPuse16, 89	tsSidebandBit8, 90
tsDigitalPuse17, 80	tsTriggered, 90
tsDigitalPuse17, 89	Unknown, 52, 53, 65, 66, 87, 90
tsDigitalPuse18, 90	unknown, 65
tsDigitalPuse19, 90	UnknownDest, 57
tsDigitalPuse2, 89	UnknownSpeed, 70
tsDigitalPuse20, 90	Unlock, 76
tsDigitalPuse21, 90	Unsigned_16bit, 61
tsDigitalPuse22, 90	Unsigned_24bit, 61
tsDigitalPuse23, 90	Unsigned_32bit, 61
tsDigitalPuse24, 90	UpdateTrigger, 85
tsDigitalPuse25, 90	USB, 53
tsDigitalPuse26, 90	USB_TARGET1, 57
tsDigitalPuse27, 90	USB_TARGET2, 57
tsDigitalPuse28, 90	USB_TARGET3, 57
tsDigitalPuse29, 90	UsbTest, 79
tsDigitalPuse3, 89	UsbVendorldEnumNet, 90
tsDigitalPuse30, 90	User_ADC_0, 68
tsDigitalPuse31, 90	User_ADC_1, 68
tsDigitalPuse4, 89	User_ADC_2, 68
tsDigitalPuse5, 89	User_ADC_3, 68
tsDigitalPuse6, 89	User_ADC_4, 68
tsDigitalPuse7, 89	User_DAC_0, 69
tsDigitalPuse8, 89	User_DAC_1, 69
tsDigitalPuse9, 89	User_DAC_2, 69
tsFeedback1, 88	UssingChamber, 67
tsFeedback10, 89	UssingClampModeEnumNet, 90
tsFeedback11, 89	UssingRail, 67
tsFeedback12, 89	UssingUnitEnumNet, 91
tsFeedback13, 89	VendorldEnumNet, 91
tsFeedback14, 89	Volt, 53, 91
tsFeedback15, 89	Voltage_3V3, 68
tsFeedback16, 89	Voltage_5V0, 68
tsFeedback17, 89	VoltageClamp, 90
tsFeedback18, 89	W2100, 79
tsFeedback19, 89	W2100_Accel_Gyro_Select_EnumNet, 91
tsFeedback2, 88	W2100DacqGroupChannelEnumNet, 91
tsFeedback20, 89	W2100DigitalSourceEnumNet, 93
tsFeedback21, 89	W2100IFB2, 67
tsFeedback22, 89	W2100Interfaceboard, 67
tsFeedback23, 89	W2100WirelessReceiver, 67, 68

W2100WirelessReceiverAnalog, 67, 68	WirelessHeadStageReservedCRE1HS1, 93
Warner, 87	WirelessHeadStageReservedCRE1HS2, 93
Warner_TEER_Machine, 79	WirelessHeadStageReservedCRE1HS3, 93
Warner_Ussing, 79	WirelessHeadStageReservedCRE1HS4, 93
WARNER_USSING_DEVICE, 62	WirelessHeadStageReservedCRE2HS1, 93
Warner_Valve_Control, 79	WirelessHeadStageReservedCRE2HS2, 93
WARNER_VALVE_CONTROL_DEVICE, 62	WirelessHeadStageReservedCRE2HS3, 93
Whole_Cell_Patch, 78	WirelessHeadStageReservedCRE2HS4, 93
WholeCellPatch, 82	WirelessHeadStageStatusRE1HS1, 92
WholeCellPatchHeadstage, 67	WirelessHeadStageStatusRE1HS2, 92
WirelessHeadStageAccDataRE1HS1, 92	WirelessHeadStageStatusRE1HS3, 92
WirelessHeadStageAccDataRE1HS2, 92	WirelessHeadStageStatusRE1HS4, 92
WirelessHeadStageAccDataRE1HS3, 92	WirelessHeadStageStatusRE2HS1, 92
WirelessHeadStageAccDataRE1HS4, 92	WirelessHeadStageStatusRE2HS2, 92
WirelessHeadStageAccDataRE2HS1, 92	WirelessHeadStageStatusRE2HS3, 92
WirelessHeadStageAccDataRE2HS2, 92	WirelessHeadStageStatusRE2HS4, 92
WirelessHeadStageAccDataRE2HS3, 92	WirelessTestAdapter, 52
WirelessHeadStageAccDataRE2HS4, 92	Work, 94
WirelessHeadStageAnalogRE1HS1, 92	WPA16, 79
WirelessHeadStageAnalogRE1HS2, 92	WPA32, 79
WirelessHeadStageAnalogRE1HS3, 92	WPA4, 79
WirelessHeadStageAnalogRE1HS4, 92	WPA8, 79
WirelessHeadStageAnalogRE2HS1, 92	WvcDisplayModeEnumNet, 94
WirelessHeadStageAnalogRE2HS2, 92	WvcValveModeEnumNet, 94
WirelessHeadStageAnalogRE2HS3, 92	Zero, 57, 71, 73, 74, 83, 86, 93
WirelessHeadStageAnalogRE2HS4, 92	MCS_ANY_BUS
WirelessHeadStageGyroDataRE1HS1, 92	Mcs::Usb, 70
WirelessHeadStageGyroDataRE1HS2, 92	MCS_CHANNELTEST_DEVICE
WirelessHeadStageGyroDataRE1HS3, 92	Mcs::Usb, 62
WirelessHeadStageGyroDataRE1HS4, 92	MCS DEVICE ANY
WirelessHeadStageGyroDataRE2HS1, 92	Mcs::Usb, 62
WirelessHeadStageGyroDataRE2HS2, 92	MCS_DEVICE_USB
WirelessHeadStageGyroDataRE2HS3, 92	Mcs::Usb, 62
WirelessHeadStageGyroDataRE2HS4, 92	MCS_DEVICE_USB_CYPRESS
WirelessHeadStageOptoStimCurrentRE1HS1, 92	Mcs::Usb, 63
WirelessHeadStageOptoStimCurrentRE1HS2, 92	MCS_ENCAPSULATOR_DEVICE
WirelessHeadStageOptoStimCurrentRE1HS3, 92	Mcs::Usb, 62
WirelessHeadStageOptoStimCurrentRE1HS4, 92	MCS_EXTERN_BC_TESTER_DEVICE
WirelessHeadStageOptoStimCurrentRE2HS1, 92	Mcs::Usb, 62
WirelessHeadStageOptoStimCurrentRE2HS2, 92	MCS_EXTERN_D_TESTER_DEVICE
WirelessHeadStageOptoStimCurrentRE2HS3, 92	Mcs::Usb, 62
WirelessHeadStageOptoStimCurrentRE2HS4, 93	MCS_FCX_DEVICE
WirelessHeadStageReservedARE1HS1, 92	Mcs::Usb, 62
WirelessHeadStageReservedARE1HS2, 92	MCS_FYI_DEVICE
WirelessHeadStageReservedARE1HS3, 92	Mcs::Usb, 62
WirelessHeadStageReservedARE1HS4, 92	MCS_GENERIC_DEVELOPMENT_DEVICE
WirelessHeadStageReservedARE2HS1, 92	Mcs::Usb, 62
WirelessHeadStageReservedARE2HS2, 92	MCS_HICLAMP_DEVICE
WirelessHeadStageReservedARE2HS3, 93	Mcs::Usb, 62
WirelessHeadStageReservedARE2HS4, 93	MCS_HLA_DEVICE
WirelessHeadStageReservedBRE1HS1, 93	Mcs::Usb, 62
WirelessHeadStageReservedBRE1HS2, 93	
	MCS_MBC08_DEVICE
WirelessHeadStageReservedBRE1HS3, 93	Mcs::Usb, 62
WirelessHeadStageReservedBRE1HS4, 93	MCS_MC_STIMULUS_DEVICE
WirelessHeadStageReservedBRE2HS1, 93	Mcs.:Usb, 62
WirelessHeadStageReservedBRE2HS2, 93	MCS_MCCARD_DEVICE
WirelessHeadStageReservedBRE2HS3, 93	Mcs::Usb, 62
WirelessHeadStageReservedBRE2HS4, 93	MCS_MEA_CLEAN_DEVICE

Mcs::Usb, 62	Mcs::Usb, 62
MCS_MEA_COAT_DEVICE	MCS_TCX_DEVICE
Mcs::Usb, 62	Mcs::Usb, 62
MCS_MEA_DEVICE	MCS_TERSENS_DEVICE
Mcs::Usb, 62	Mcs::Usb, 62
MCS_MEA_IMPEDANCE_DEVICE	MCS UNDEFINED BUS
Mcs::Usb, 62	Mcs::Usb, 70
MCS_MEA_SWITCH_DEVICE	MCS_USB_BUS
Mcs::Usb, 62	Mcs::Usb, 70
MCS MEASURETABLE DEVICE	McsBus
-	
Mcs::Usb, 62	CPPC DeviceNet, 424
MCS_MEAUSB_DEVICE	CPPS_DeviceNet, 433
Mcs::Usb, 62	CRoboDeviceNet, 486
MCS_NF_GEN_DEVICE	MCSBUS0
Mcs::Usb, 62	Mcs::Usb, 54
MCS_OCTOPOT_DEVICE	MCSBUS1
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_OKUVISION_STIMULATOR_DEVICE	Mcs::Usb, 53
Mcs::Usb, 62	MCSBUS10
MCS_PATCHSERVER_DEVICE	FirmwareDestinationNames, 692
Mcs::Usb, 62	Mcs::Usb, 54
MCS_PATHIDENT_DEVICE	MCSBUS11
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_PCI_BUS	Mcs::Usb, 54
Mcs::Usb, 70	MCSBUS12
MCS_PCX_DEVICE	FirmwareDestinationNames, 692
Mcs::Usb, 62	Mcs::Usb, 54
MCS_PEDOTER_DEVICE	MCSBUS13
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_PERISTALTIC_PUMP_DEVICE	Mcs::Usb, 54
Mcs::Usb, 62	MCSBUS14
MCS_PGA_DEVICE	Mcs::Usb, 54
Mcs::Usb, 62	MCSBUS15
MCS_PPC_DEVICE	Mcs::Usb, 54
Mcs::Usb, 62	MCSBUS2
MCS_PPS5_DEVICE	FirmwareDestinationNames, 692
Mcs::Usb, 62	Mcs::Usb, 53
MCS_PPS_DEVICE	MCSBUS3
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_RETINA_AMS_DONGLE	Mcs::Usb, 53
Mcs::Usb, 62	MCSBUS4
MCS_RETINA_LED_DEVICE	FirmwareDestinationNames, 692
Mcs::Usb, 62	Mcs::Usb, 53
MCS_ROBO_DEVICE	MCSBUS5
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_ROBOINJECT_DEVICE	Mcs::Usb, 53
Mcs::Usb, 62	MCSBUS6
MCS_ROBOOCYTE2_DEVICE	FirmwareDestinationNames, 692
Mcs::Usb, 62	Mcs::Usb, 54
MCS_SAFEIS_DEVICE	MCSBUS7
Mcs::Usb, 62	FirmwareDestinationNames, 692
MCS_SMARTIMPLANT_DEVICE	Mcs::Usb, 54
Mcs::Usb, 62	MCSBUS8
MCS_SOFTWARE_DONGLE_DEVICE	FirmwareDestinationNames, 693
Mcs::Usb, 62	Mcs::Usb, 54
MCS_STG_DEVICE	MCSBUS9
Mcs::Usb, 62	FirmwareDestinationNames, 693
MCS_SW2TO64_DEVICE	Mcs::Usb, 54
INIOO_ONNETOUT_DENIOE	IVICOCOD, 04

McsBus_MotorControl	Mcs::Usb, 78
CPeristalticPumpDeviceNet, 409	MEA1060
CPPCDeviceNet, 424	Mcs::Usb, 78
CPPS_DeviceNet, 433	MEA120
CRoboDeviceNet, 486	Mcs::Usb, 52
CRoboFluidDeviceNet, 490	MEA2100
McsBus_Sensor	Mcs::Usb, 78
CPPCDeviceNet, 424	Mea2100
CPPS_DeviceNet, 433	Mcs::Usb, 67
McsBus_VoltageMode	MEA2100 256
CFluidControlDeviceNet, 149	Mcs::Usb, 78
McsBus XY	Mea2100 256
CRoboDeviceNet, 483	Mcs::Usb, 68
McsBus ZI	MEA2100_256DacqGroupChannelEnumNet
CRoboDeviceNet, 483	Mcs::Usb, 71
McsBusTypeEnumNet	MEA2100_256DigitalSourceEnumNet
Mcs::Usb, 70	Mcs::Usb, 71
McsUsbDeviceStateEvent	Mea2100_256Headstage
CMcsUsbDeviceStatePushFunctionNet, 298	Mcs::Usb, 67
CMcsUsbDeviceStatePushNet, 299	Mea2100 256Interfaceboard
McsUsbSpeedEnumNet	_
Mcs::Usb, 70	Mcs::Usb, 67 MEA2100 32
MCU1	-
	Mcs::Usb, 78
FirmwareDestinationNames, 693	MEA2100_Lite
Mcs::Usb, 53	Mcs::Usb, 78
ME128	Mea2100_Lite
Mcs::Usb, 78	Mcs::Usb, 68
ME16	MEA2100_Mini
Mcs::Usb, 78	Mcs::Usb, 78
ME2100	MEA2100_Mini_Usb_develop
Mcs::Usb, 78	Mcs::Usb, 78
Me2100_32PICiCE40	MEA2100BetaScreen
Mcs::Usb, 81	Mcs::Usb, 78
Me2100_32PICiCE40Headstage	Mea2100BetaScreen
Mcs::Usb, 67	Mcs::Usb, 81
Me2100_32Xilinx	Mea2100BetaScreenHeadstage
Mcs::Usb, 81	Mcs::Usb, 67
Me2100_32XilinxHeadstage	MEA2100DigitalSourceEnumNet
Mcs::Usb, 67	Mcs::Usb, 73
Me2100Graphene16_32	Mea2100Headstage
Mcs::Usb, 82	Mcs::Usb, 66
Me2100Graphene16_32Headstage	Mea2100Interfaceboard
Mcs::Usb, 67	Mcs::Usb, 66
Me2100Interfaceboard	Mea2100LiteHeadstage
Mcs::Usb, 67	Mcs::Usb, 67
Me2100InvitroSignalCollectorUnit	Mea2100Mini120
Mcs::Usb, 67	Mcs::Usb, 82
Me2100InvivoSignalCollectorUnit	Mea2100Mini120Headstage
Mcs::Usb, 67	Mcs::Usb, 67
Me2100UPA32	Mea2100Mini60ECP5
Mcs::Usb, 81	Mcs::Usb, 82
Me2100UPA32Headstage	Mea2100Mini60ECP5Headstage
Mcs::Usb, 67	Mcs::Usb, 67
ME256	Mea2100Mini60PICiCE40
Mcs::Usb, 78	Mcs::Usb, 81
ME32	Mea2100Mini60PICiCE40Headstage
Mcs::Usb, 78	Mcs::Usb, 67
	Mea2100MultiwellIFB2
ME64	IVICAL I UUIVIUILIWEIIIFDL

Mcs::Usb, 67	Mcs::Usb, 60
Mea2100STG	mkfilter
Mcs::Usb, 66	mkfilterNet, 699
MEA252	mkfilter_coef_in_one_set
Mcs::Usb, 52	mkfilterNet, 699
MEA256	mkfilter_highpass_coeff
Mcs::Usb, 78	mkfilterNet, 700
MEA2x32	mkfilter_highpass_frequency_from_coeff
Mcs::Usb, 52	mkfilterNet, 700
MEA2x60	mkfilter_highpass_frequency_from_k
Mcs::Usb, 52	mkfilterNet, 700
MEA32	mkfilter highpass k
Mcs::Usb, 52	mkfilterNet, 700
MEA60	mkfilter_MCS
Mcs::Usb, 52	mkfilterNet, 700
MEA_2_252_2	mkfilter_MCS_k
Mcs::Usb, 52	mkfilterNet, 700, 701
MEA_2_252_2_6Well	mkfilter normalize coeffs int
Mcs::Usb, 52	mkfilterNet, 701
MEA_2_252_2_9Well	mkfilter normalize coeffs short
Mcs::Usb, 52	mkfilterNet, 701
MEA_2_252_2_Test	mkfilter_normalize_scale_coeffs_int
Mcs::Usb, 52	mkfilterNet, 701
MEA Clean	mkfilter_scale_coef_in_one_set
Mcs::Usb, 78	mkfilterNet, 701
MEA_Coat	mkfilterNet, 699
Mcs::Usb, 78	mkfilter, 699
MEA_Impedance	mkfilter_coef_in_one_set, 699
Mcs::Usb, 78	mkfilter_highpass_coeff, 700
MEA_Sanofi	mkfilter_highpass_frequency_from_coeff, 700
Mcs::Usb, 78	mkfilter_highpass_frequency_from_k, 700
MEA_Switch	mkfilter_highpass_k, 700
Mcs::Usb, 78	mkfilter_MCS, 700
MEA_Switch_2_1	mkfilter_MCS_k, 700, 701
Mcs::Usb, 78	mkfilter_normalize_coeffs_int, 701
MEA_Switch_4_2	mkfilter_normalize_coeffs_short, 701
Mcs::Usb, 78	mkfilter_normalize_scale_coeffs_int, 701
MeaAudioFunctionNet	mkfilter_scale_coef_in_one_set, 701
CMeaDeviceNet, 361	mIMEA60
MeaDigitalDataFunctionNet	Mcs::Usb, 74
CMeaDeviceNet, 361	mlUnknown
MeaFeedbackFunctionNet	Mcs::Usb, 74
CMeaDeviceNet, 361	ModulA_ADC0
MeaLayoutEnumNet	Mcs::Usb, 68
Mcs::Usb, 74	ModulA_ADC1
Measure	Mcs::Usb, 68
CPathIdentDeviceNet, 406	ModulA_ADC2
MeasureReservoir	Mcs::Usb, 68
CPPCFunctionNet, 431	ModulA_ADC3
MeasureTable	Mcs::Usb, 68
Mcs::Usb, 79	ModulA_DAC0
MeasuringOnly	Mcs::Usb, 69
HeadStageIDType, 694	ModulA_DAC1
MeFunctionNet	Mcs::Usb, 69
CMeaDeviceNet, 361	ModulB_ADC0
MicroAmpere	Mcs::Usb, 68
Mcs::Usb, 60	ModulB_ADC1
MilliDegreeCelsius	Mcs::Usb, 68

ModulB_ADC2	Multiwell96
Mcs::Usb, 68	Mcs::Usb, 52
ModulB_ADC3	Multiwell_ICC
Mcs::Usb, 68	Mcs::Usb, 78
ModulB_DAC0	Multiwell_MEA_Mini
Mcs::Usb, 69	Mcs::Usb, 78
ModulB_DAC1	MultiwellDigitalSourceEnumNet
Mcs::Usb, 69	Mcs::Usb, 74
ModulC_ADC0	MultiwellHeadstage
Mcs::Usb, 68	Mcs::Usb, 66
ModulC_ADC1	MultiwellInterfaceboard
Mcs::Usb, 68	Mcs::Usb, 66
ModulC_ADC2	MultiwellMini
Mcs::Usb, 68	Mcs::Usb, 82
ModulC_ADC3	MultiwellMiniHeadstage
Mcs::Usb, 68	Mcs::Usb, 67
ModulC DAC0	MultiwellOptoStim
Mcs::Usb, 69	Mcs::Usb, 77
ModulC DAC1	MultiwellPlateTypeEnumNet
Mcs::Usb, 69	Mcs::Usb, 75
ModulD ADC0	Mux
Mcs::Usb, 69	Mcs::Usb, 63
ModulD ADC1	MuxOtherDevice
Mcs::Usb, 69	Mcs::Usb, 63
ModulD ADC2	MwPollStatusEvent
Mcs::Usb, 69	CStg200xDownloadNet, 573
ModulD ADC3	Cotg200xBowindaditot, 070
Mcs::Usb, 69	Nanion
ModulD DAC0	Mcs::Usb, 87
Mcs::Usb, 69	NanoAmpere
ModulD DAC1	Mcs::Usb, 60
Mcs::Usb, 69	NanoVolt
Monitor	Mcs::Usb, 60
Mcs::Usb, 84	NCBathCurrent
MoveAbs	Mcs::Usb, 58
CRoboDeviceNet, 478, 479	NCCol2Current
MoveAbsl	Mcs::Usb, 58
CRoboStatorDeviceNet, 495	NChipTemperature
MoveAbsXY	Mcs::Usb, 58
	Neptun
CRoboStatorDeviceNet, 495	
	Mcs::Usb. 79
MoveAbsZ	Mcs::Usb, 79 NeuroChip
CRoboStatorDeviceNet, 495, 496	NeuroChip
CRoboStatorDeviceNet, 495, 496 Movement	NeuroChip Mcs::Usb, 79
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81	NeuroChip Mcs::Usb, 79 NeurochipConfig
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage CMcsUsbNet, 330	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage CMcsUsbNet, 330 MultibootSelectImage	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage CMcsUsbNet, 330 MultibootSelectImage CMcsUsbNet, 330	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65 nMos36LN
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage CMcsUsbNet, 330 MultibootSelectImage CMcsUsbNet, 330 MultibootSelectImage CMcsUsbNet, 330 MultibiootSelectImage CMcsUsbNet, 330 MultiFileMode	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65 nMos36LN Mcs::Usb, 65
CRoboStatorDeviceNet, 495, 496 Movement	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65 nMos36LN Mcs::Usb, 65 nMos36LN Mcs::Usb, 65
CRoboStatorDeviceNet, 495, 496 Movement Mcs::Usb, 81 Multiboot Mcs::Usb, 79 MultibootGetCypressImageId CMcsUsbNet, 330 MultibootGetImageId CMcsUsbNet, 330 MultibootGetSelectedImage CMcsUsbNet, 330 MultibootSelectImage CMcsUsbNet, 330 MultibiootSelectImage CMcsUsbNet, 330 MultiFileMode Mcs::Usb, 84 Multiwell	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65 nMos36LN Mcs::Usb, 65 nMos64LN Mcs::Usb, 65
CRoboStatorDeviceNet, 495, 496 Movement	NeuroChip Mcs::Usb, 79 NeurochipConfig Mcs::Usb, 79 NF_Gen Mcs::Usb, 79 NineWell Mcs::Usb, 86 nMos16LV Mcs::Usb, 65 nMos32LV Mcs::Usb, 65 nMos36LN Mcs::Usb, 65 nMos36LN Mcs::Usb, 65

Mcs::Usb, 52 NOT_CONNECTED Mcs::Usb, 55 NOT_CONNECTED Mcs::Usb, 55 NOtApplicable Mcs::Usb, 52 NoUnit Mcs::Usb, 60 NTC10K Mcs::Usb, 87 NUICommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfIchannels HeadStageIDType, 695 NumberOfStimulationChannels HeadStageIDType, 696 NumCoelfSter CoreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 79 Off Mcs::Usb, 80 OK Mcs::Usb, 80 ONGelAnalogThresholdicub Mcs::Usb, 94 OnedeAnalogThresholdicub CWarnerValveControllerDeviceNet, 660 OnGelAnalogThresholdicub CWarnerValveControllerDeviceNet, 660 OnGelAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGelAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGelDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGelDigitalOutPortVa	None	CWarnerValveControllerDeviceNet, 660
Mcs::Usb, 55 NotApplicable Mcs::Usb, 52 NotApplicable Mcs::Usb, 50 NotApplicable Mcs::Usb, 50 NotTolk Mcs::Usb, 57 NotCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStagelDType, 695 NumberOfAnalogChannels HeadStagelDType, 695 NumberOfStimulationChannels HeadStagelDType, 696 NumcoelSets CreateFilterNet, 118 CotociceGroupChannelInfoTemplateNet > , 127 Off Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 77 Off Mcs::Usb, 80 OK Mcs::Usb, 91 OK Mcs::Usb, 95 OKStgPoliStatus Mcs::Usb, 95 ONTableIentryChanged CWarnerValveControllerDeviceNet, 660 OKGetAnalogThreshold-ligh CWarnerValveControllerDeviceNet, 660 OKGetAlanalogThreshold-ligh CWarnerValveControllerDeviceNet, 660 OKGetDigitalOrDitrollerDeviceNet, 660 OKGetDigitalOrDitrollerDeviceNet, 660 OKGetDigitalOrDitrollerDeviceNet, 660 OKGetDigitalPorticute CWarnerValveControllerDeviceNet, 660 OKGetDigitalPorticuteNet CWarnerValveControllerDeviceNet, 660 OKGetDigitalPorticuteNet CWarnerValveControllerDevic	Mcs::Usb, 52, 63, 77, 85, 90, 91	OnGetPlateClampStateByHeadstage
NOT_CONNECTED Mcs::Usb, 55 NotApplicable Mcs::Usb, 52 NoUnint Mcs::Usb, 60 NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageiDType, 695 NumberOfChanelEnumTemplateNet < Darcy GroupChannelIntoTemplateNet < Darcy GroupChannelIntoTemplateNet < Darcy GroupChannelIntoTemplateNet < Darcy Mcs::Usb, 77 NumberOfStimulationChannels HeadStageiDType, 695 NumcoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 79, 81, 91 off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 91 OK OREValveDigitalInPort CWarnerValveControllerDeviceNet, 660 OnicetAnalogitalProviceNet, 660 OnicetAnalogitalProviceNet, 660 OnicetAnalogitalProviceNet, 660 OnicetAnalogitalProviceNet, 660 OnicetAnalogitalProviceNet, 660 OnicetCurrentNumberOtValves CWarnerValveControllerDeviceNet, 660 OnicetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnicetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnicetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnicetDigitalPortivalves CWarnerValveControllerDeviceNet, 660 OnicetDigitalPortivalve CWarnerValveControllerDeviceNet, 660 OnicetDi	Normal	CMultiwellCallbackFunctionNet, 383
Mcs::Usb, 55 NotApplicable Mcs::Usb, 52 NotDrit Mcs::Usb, 50 NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStagelDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet < DackgroupChannelEnumTemplateNet > 127 NumberOfStimulationChannels HeadStagelDType, 695 NumberOfStimulationChannels HeadStagelDType, 696 NumCoerSets CCreateFitterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 ONGenameIData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetGDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetGDigi	Mcs::Usb, 75	OnGetTableNamebyIndex
NotApplicable Mcs::Usb, 52 NoUnit Mcs::Usb, 60 NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfChannels CDeviceGroupChannelInumTemplateNet < package of the provided and th	NOT_CONNECTED	CWarnerValveControllerDeviceNet, 660
Mcs::Usb, 52 NoUnit Mcs::Usb, 60 NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnanels HeadStageIDType, 695 NumberOfShannels HeadStageIDType, 695 NumberOfShannelInfoTemplateNet < GroupChannelEnumTemplateNet > 127 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 79 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 56 Okuvision_Stimulator Mcs::Usb, 78 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnGeltAcitiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetDigitalOrtortlaye CWarnerValveControllerDeviceNet, 660 OnGetOltigitalOrtortlaye	Mcs::Usb, 65	OnGetValveActive
Mcs::Usb, 52 NoUnit Mcs::Usb, 60 NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfStimulationChannels HeadStageIDType, 695 NumberOfStimulationChannels HeadStageIDType, 696 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 56 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 78 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalveControllerDeviceNet, 660 OnGetAlactiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLing CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogUtlage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetGloigitalOurPortValve CWarnerValveControllerDeviceNet, 660 OnGetOptiglatOurPortValve CWarnerValveControllerDeviceNet, 660 OnGet	NotApplicable	CWarnerValveControllerDeviceNet, 661
Mcs::Usb, 60 NTG10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStagelDType, 695 NumberOfAnanels CDeviceGroupChannelInfoTemplateNet < DackgroupChannelEnumTemplateNet > 127 NumberOfStimulationChannels HeadStagelDType, 696 NumCoefSets CCreateFilterNet, 118 Cctopot Mcs::Usb, 77 Off Ctotopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAnatogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogVitage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVitage CWarnerValveControllerDeviceNet, 660 OnGetGligitalOutPortValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNe	Mcs::Usb, 52	OnGetValveBoardRevision
Mcs::Usb, 60 NTG10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStagelDType, 695 NumberOfAnanels CDeviceGroupChannelInfoTemplateNet < DackgroupChannelEnumTemplateNet > 127 NumberOfStimulationChannels HeadStagelDType, 696 NumCoefSets CCreateFilterNet, 118 Cctopot Mcs::Usb, 77 Off Ctotopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAnatogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogVitage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVitage CWarnerValveControllerDeviceNet, 660 OnGetGligitalOutPortValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNe	NoUnit	CWarnerValveControllerDeviceNet, 661
NTC10K Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet <	Mcs::Usb, 60	
Mcs::Usb, 87 NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet < DackgroupChannelEnumTemplateNet >, 127 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAlaveNorder OneGetAlaveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetCluveNetControllerDeviceNet, 660 OnGetCluveNetControllerD		•
NullCommand CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet <		
CRoboDeviceNet, 480 NumberOfAnalogChannels HeadStageIDType, 695 NumberOfChannels CDeviceGroupChannelInumTemplateNet < Dacq- GroupChannelEnumTemplateNet > , 127 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControll		
NumberOfAnalogChannels HeadStagelDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet < Dacq- GroupChannelEndTemplateNet > , 127 NumberOfStimulationChannels HeadStagelDType, 696 NumCoefSets CCreateFilterNet, 118 Cotopot Mcs::Usb, 77 Off Mcs::Usb, 77 Off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 85 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 78 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve		
HeadStagelDType, 695 NumberOfChannels CDeviceGroupChannelInfoTemplateNet < Dacq- GroupChannelEnumTemplateNet >, 127 NumberOfStimulationChannels HeadStagelDType, 696 NumCoefSets CCreateFilterNet, 118 Cctopot Mcs::Usb, 77 Off Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 85 OnChannelData Mcs::Usb, 86 OKuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdI-ligh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdI-low CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdI-low CWarnerValveControllerDeviceNet, 660 OnGetAnalogTortollerDeviceNet, 660 OnGetAnalogTortollerDeviceNet, 660 OnGetAnalogTortollerDeviceNet, 660 OnGetDigitalOutPortNalve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortNalve CWarnerVal		•
NumberOfChannels CDeviceGroupChannelInfoTemplateNet < Dacq- GroupChannelEnumTemplateNet > 127 NumberOfStimulationChannels HeadStagelDType, 696 NumCoefSets CCreateFilterNet, 118 Cctopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogValveControllerDeviceNet, 660 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortIverted CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalIntorer CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalMode CWarnerValveControllerDeviceNet, 662 OnMsulsbDeviceState Mcs::Usb, 94 OnMcs::Usb, 95 OnMsulsbDeviceState Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnDupdateFirmwareProgress Mcs::Usb, 95 OnpenPipe CogenClamp Mcs::Usb, 96 OpenClamp Mcs::Usb, 96 OpenPipe CGenericDevelopDeviceNet, 162 OpenPipe CGenericDevelopDeviceNet, 389 Operator=		
CDeviceGroupChannelEnumTemplateNet >, 127 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 85 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 OnGetAalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDevi	- · · ·	
GroupChannelEnumTemplateNet >, 127 NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 off Mcs::Usb, 80 OK Mcs::Usb, 85 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveContr		
NumberOfStimulationChannels HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Cctopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnEtror Mcs::Usb, 94 One Mcs::Usb, 95 OnEtror Mcs::Usb, 94 OnGetAnalogThresholdLidgh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 One CWarnerValveControllerDeviceNet, 660 One CWarnerValveControllerDeviceNet, 660 One CWarnerValveControl	·	
HeadStageIDType, 696 NumCoefSets CCreateFilterNet, 118 Cotopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLiow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection C	·	•
NumCoefSets		•
CCreateFilterNet, 118 CSCUFunctionNet, 515 OnlsValveDigitalInInverted CWarnerValveControllerDeviceNet, 662 OnlsValveOpen CWarnerValveControllerDeviceNet, 662 OnlsValveOpen CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInAnalogMode CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalMode CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalMode CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalMode CWarnerValveControllerDeviceNet, 662 OnlsValveOpenInDigitalMode CWarnerValveControllerDeviceNet, 662 OnloopticeArrivalRemoval Mcs::Usb, 65 OnCettopic Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 Onstg200xErarlandler Mcs::Usb, 95 Onstg20loxErarlandler Mcs::Usb, 95 Onstg20loxErarlandler Mcs::Usb, 95 OnstgPollStatus Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CWarnerValveControllerDeviceNet, 660 OnGetEdigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OngetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OngetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OngetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OngetDigitalOutPortValve OpenPipe CGenericDevelopDeviceNet, 162 OpenPipe	5	
Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValve		
Octopot Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 Onerror Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection	Corealermen, 116	
Mcs::Usb, 77 Off Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnEtror Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection CWarnerValveControllerDeviceNet, 660	Octopot	3
Off	•	
Mcs::Usb, 79, 81, 91 Off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnsetDigitalPortDirection CWarnerValveContr		•
off Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 O		
Mcs::Usb, 80 OK Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetAralogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection OnGetDigitalPortDirection OnGetDigitalPortDirection OnGetDigitalPortDirection OnGetDigital		
OK Mcs::Usb, 65 Ckuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnJedateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 96 Open Mcs::Usb, 96 OpenClamp Mcs::Usb, 76 OpenClamp Mcs::Usb, 96 OpenClamp Mcs::Usb, 90 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 operator=		CWarnerValveControllerDeviceNet, 662
Mcs::Usb, 65 Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStgPollStatus OnJupdateFirmwareProgress Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660		OnIsValveOpenInDigitalMode
Okuvision_Stimulator Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalo		CWarnerValveControllerDeviceNet, 662
Mcs::Usb, 78 OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnUpdateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 95 Open Mcs::Usb, 95 OpenClamp Mcs::Usb, 96 OpenClamp Mcs::Usb, 90 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 operator=		OnMcsUsbDeviceState
OnChannelData Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDi	Okuvision_Stimulator	Mcs::Usb, 94
Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection OndetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OndetDigitalPortDirection	•	OnMcsUsbDeviceStateCallback
Mcs::Usb, 94 OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPo	OnChannelData	Mcs::Usb, 95
OnDeviceArrivalRemoval Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnStg200xDataHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnStgPollStatus Mcs::Usb, 95 OnStgPollStatus Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection Onget	Mcs::Usb, 94	
Mcs::Usb, 94 One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=	OnDeviceArrivalRemoval	
One Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPlateClamp CMultiwellDeviceNet, 389 OnStg200xErrorHandler Mcs::Usb, 95 OnStg200xErrorHandler Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 662 OnUpdateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 95 OnUpdateFirmwareProgress Mcs::Usb, 95 Open Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::	Mcs::Usb, 94	
Mcs::Usb, 57, 71, 73, 74, 83, 86, 93 OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnStg200xErrorHandler Mcs::Usb, 95 OnTableEntryChanged CWarnerValveControllerDeviceNet, 662 OnUpdateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 95 OnUpdateFirmwareProgress OnEstiveControllerDeviceNet, 660 Open Mcs::Usb, 95 OnUpdateFirmwareProgress OnUpdateFirmwareProgress OnUpdateFirmwareProgress OnUpdateFirmwareProgress OnUpdateFirmwareProgress OnUpdateFirmwareProgress OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 96 Open Mcs:	One	
OnError Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=	Mcs::Usb, 57, 71, 73, 74, 83, 86, 93	
Mcs::Usb, 94 OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigita		
OnGetActiveRunningTableNumber CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection OnGetDigitalPortDirection		
CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnTableEntryChanged CWarnerValveControllerDeviceNet, 660 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Mcs::Usb, 95 OpenClamp Mcs::Usb, 76 OpenClamp Mcs::Usb, 90 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		
OnGetAnalogThresholdHigh CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=	•	
CWarnerValveControllerDeviceNet, 660 OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OnUpdateFirmwareProgress Mcs::Usb, 95 OnUpdateFirmwareStatusChange Mcs::Usb, 95 Open Copen Mcs::Usb, 76 OpenClamp Mcs::Usb, 90 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		
OnGetAnalogThresholdLow CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		
CWarnerValveControllerDeviceNet, 660 OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		
OnGetAnalogVoltage CWarnerValveControllerDeviceNet, 660 OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		
CWarnerValveControllerDeviceNet, 660 Open OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OpenClamp Mcs::Usb, 76 OpenClamp Mcs::Usb, 90 OpenPipe CGenericDevelopDeviceNet, 162 CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		OnUpdateFirmwareStatusChange
OnGetAvailableHeadstages CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=	• •	Mcs::Usb, 95
CSCUFunctionNet, 515 OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		Open
OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPlateClamp CMultiwellDeviceNet, 389 Operator=		Mcs::Usb, 76
OnGetCurrentNumberOfValves CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 operator=		OpenClamp
CWarnerValveControllerDeviceNet, 660 OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPipe CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 operator=		·
OnGetDigitalOutPortValve CWarnerValveControllerDeviceNet, 660 OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 CGenericDevelopDeviceNet, 162 OpenPlateClamp CMultiwellDeviceNet, 389 operator=	CWarnerValveControllerDeviceNet, 660	
CWarnerValveControllerDeviceNet, 660 OpenPlateClamp OnGetDigitalPortDirection CWarnerValveControllerDeviceNet, 660 CWarnerValveControllerDeviceNet, 660 OpenPlateClamp CMultiwellDeviceNet, 389 operator=	OnGetDigitalOutPortValve	•
OnGetDigitalPortDirection	_	•
CWarnerValveControllerDeviceNet, 660 operator=		·
oporator—		
	OnGetDisplayMode	ορθιαίοι=

DeviceIdNet, 681	PicoAmpere
OpticalStimulation	Mcs::Usb, 60
HeadStageIDType, 694	Plate_24W030MGA
Order	Mcs::Usb, 75
CCreateFilterNet, 118	Plate_24W300_30_1152GBA
CFilterPropertyNet, 142	Mcs::Usb, 75
Output	Plate_24W300_30GBA
Mcs::Usb, 76	Mcs::Usb, 75
De election of October 40 and the	Plate_24W300_30GBB
PacketFrameContextGroup	Mcs::Usb, 75
Mcs::Usb, 61, 71, 82, 93	Plate_24W300_30GMA
PatchServAdcModeEnumNet	Mcs::Usb, 75
Mcs::Usb, 75	Plate_24W700_100FMA
PatchServer	Mcs::Usb, 75
Mcs::Usb, 79	Plate_24W700_100FMB
PathIdent	Mcs::Usb, 75
Mcs::Usb, 79	Plate_24W700_100FMC
PatternListStart	Mcs::Usb, 75
COctoPotDeviceNet, 399	Plate_24W700_100PBA
PC	Mcs::Usb, 75
Mcs::Usb, 94	Plate_72W500_100FMA
PCI	Mcs::Usb, 75
Mcs::Usb, 91	Plate_72W500_100PMA
PCX	Mcs::Usb, 75
Mcs::Usb, 77	Plate_96W300_80_1152FMA
PeriodicPulse	Mcs::Usb, 75
Mcs::Usb, 63	Plate_96W400_80_1152FMB
PeristalticPump	Mcs::Usb, 75
Mcs::Usb, 78	Plate_96W700_100FMA
PGA	Mcs::Usb, 75
Mcs::Usb, 77	Plate_96W700_100FMB
PIC	Mcs::Usb, 75
FirmwareDestinationNames, 693	Plate_96W700_100GBA
Mcs::Usb, 55 PIC10	Mcs::Usb, 75
	Plate_96W700_100GBB
Mcs::Usb, 55 PIC11	Mcs::Usb, 75
	Plate_96W700_100GBC
Mcs::Usb, 55	Mcs::Usb, 75
PIC12	Plate_96W700_100GBD
Mcs::Usb, 55	Mcs::Usb, 75
PIC2	Plate_96W700_100GMA
FirmwareDestinationNames, 693	Mcs::Usb, 75
Mcs::Usb, 55	Plate_Dummy
PIC3	Mcs::Usb, 75
FirmwareDestinationNames, 693	Plate_Dummy_126
Mcs::Usb, 55	Mcs::Usb, 75
PIC4	PlateClampEnumNet
FirmwareDestinationNames, 693	Mcs::Usb, 76
Mcs::Usb, 55	PlateClampLockEnumNet
PIC5	Mcs::Usb, 76
Mcs::Usb, 55	PlusMinus10Volts
PIC6	Mcs::Usb, 53
Mcs::Usb, 55	PlusMinus2Comma5Volts
PIC7	Mcs::Usb, 53
Mcs::Usb, 55	PlusMinus5Volts
PIC8	Mcs::Usb, 53
Mcs::Usb, 55 PIC9	PollStatusEvent
	CStimulusFunctionNet, 584
Mcs::Usb, 55	

CW2100_StimulatorFunctionNet, 628	Product
PortDirectionEnumNet	CMcsUsbListEntryNet, 314
Mcs::Usb, 76	ProductIdEnumNet
Pos900	Mcs::Usb, 76
Mcs::Usb, 79	Program
PositionBase	CProgramPressureCurveNet, 440
Mcs::Usb, 79	PT100
PositionIIBase	Mcs::Usb, 87
Mcs::Usb, 79	PT1000
PositionIICentralUnit	Mcs::Usb, 87
Mcs::Usb, 79	PulseGenerator
PositionImp	CW2100 FunctionNet, 621
Mcs::Usb, 79	Mcs::Usb, 57, 71, 74, 75, 83, 86, 93
PostCommaA	PulseGenerator_Mode_EnumNet
CFilterCoefficientsNet::s_FilterAttributesNet, 710	Mcs::Usb, 79
Mcs::Usb, 65	PumpOff
PostCommaB	CRoboFluidDeviceNet, 489
CFilterCoefficientsNet::s FilterAttributesNet, 710	PumpOn
Mcs::Usb, 65	CRoboFluidDeviceNet, 489
PowerChip	Chobol luidbevicenet, 489
CCMOSMea_FunctionNet, 111	QueryTriggerstatus
PowerHS	CStg200xDownloadNet, 571
	Gotg200x20Willoadivot, 071
CSCUFunctionNet, 515	RampStart
PP_Pump_Mode_Type_EnumNet	COctoPotDeviceNet, 399
Mcs::Usb, 76	rawdata
PPC	Mcs::Usb, 86
Mcs::Usb, 78	RC
PPCFunction CREAT AND A 10 C	Mcs::Usb, 66
CPPCDeviceNet, 424	rc100mAh
PPS2	Mcs::Usb, 70
Mcs::Usb, 78	rc200mAh
PPS4plus1	Mcs::Usb, 70
Mcs::Usb, 78	rc300mAh
PPS5	Mcs::Usb, 70
Mcs::Usb, 78	rc30mAh
PPS5_DIG	Mcs::Usb, 70
Mcs::Usb, 78	rcGreater300mAh
PPS_Function	Mcs::Usb, 70
CPPS_DeviceNet, 433	Read
PreCommaA	
CFilterCoefficientsNet::s_FilterAttributesNet, 710	CExternDTesterDeviceNet, 134
Mcs::Usb, 65	Read2
PreCommaB	CExternDTesterDeviceNet, 134
CFilterCoefficientsNet::s_FilterAttributesNet, 710	ReadBlockFromFlash
Mcs::Usb, 65	CMcsUsbFactoryNet, 305
PrepareAndAppendData	ReadBlockFromIFBGlobalEEprom
CStg200xDownloadNet, 569	CMcsUsbFactoryNet, 305
CStimulusFunctionNet, 579	ReadBlockFromNVMEM
PrepareAndSendData	CMcsUsbFactoryNet, 305
CStg200xDownloadNet, 570	ReadClipping
CStimulusFunctionNet, 580	CLIH3DeviceNet, 201
PrepareChannelData	ReadEepromRegisterPreconfig
CDigOutStimulatorFunctionNet, 130	CMcsUsbNet, 330, 331
PrepareData	ReadPipe
CStimulusFunctionNet, 581	CGenericDevelopDeviceNet, 162
CW2100_StimulatorFunctionNet, 626	ReadRegister
PrepareDataSync	CMcsUsbNet, 331
CW2100_StimulatorFunctionNet, 626	ReadRegister32
	CMcsUsbNet, 331

B 18 1 2 TI 81 2	
ReadRegisterTimeSlot	Restart
CMcsUsbNet, 331	Mcs::Usb, 80
ReadUARTData	Retina_LED
CLIH3DeviceNet, 201	Mcs::Usb, 77
Receive	RetriggerActionEnumNet
CSerialPortNet, 520	Mcs::Usb, 80
ReceiveString	RFFunction
CSerialPortNet, 520	CPositionIIDeviceNet, 420
Rectangle	Rising
Mcs::Usb, 88	Mcs::Usb, 64
Ref16	RoboCurrentModeEnumNet
Mcs::Usb, 80	Mcs::Usb, 80
Ref24	RoboDacq
Mcs::Usb, 80	CHiClampDeviceNet, 184
Ref32	RoboDevice
Mcs::Usb, 80	CSafeISDeviceNet, 501
Ref8	
	RoboError_AnotherMaster
Mcs::Usb, 80	CRoboDeviceNet, 483
Reference	RoboError_Base
Mcs::Usb, 81	CRoboDeviceNet, 484
ReferenceElectrodeModeEnumNet	RoboError_CannotEscapeEndSwitch
Mcs::Usb, 80	CRoboDeviceNet, 484
ReferenceElectrodeSwitchPositionEnumNet	RoboError_CommandAlreadyInProgress
Mcs::Usb, 80	CRoboDeviceNet, 484
RegisterHigh	RoboError_CommandNotPossible
Mcs::Usb, 63	CRoboDeviceNet, 484
RegisterLow	RoboError_CommunicationTimeout
Mcs::Usb, 63	CRoboDeviceNet, 484
Regular	RoboError_DacqNotReady
Mcs::Usb, 87	CRoboDeviceNet, 484
RemoveSoftwareKey	RoboError_DLLMovementTimeout
CMcsUsbNet, 331	CRoboDeviceNet, 484
Renesas	RoboError FindReferenceMethod
Mcs::Usb, 90	CRoboDeviceNet, 484
RescanHeadstage	RoboError GilsonCommandPending
CMcsUsbNet, 331	CRoboDeviceNet, 484
Reserved1	RoboError GilsonTimeout
Mcs::Usb, 87	-
•	CRoboDeviceNet, 484
Reserved2	RoboError_GilsonWrondID
Mcs::Usb, 87	CRoboDeviceNet, 485
Reserved3	RoboError_McsBus_UnknownCommand
Mcs::Usb, 87	CRoboDeviceNet, 485
Reserved4	RoboError_NoEndSwitch
Mcs::Usb, 87	CRoboDeviceNet, 485
Reserved5	RoboError_NoMoreData
Mcs::Usb, 87	CRoboDeviceNet, 485
ResetAdcOffset	RoboError_NoReference
COctoPotDeviceNet, 399	CRoboDeviceNet, 485
ResetChannelmap	RoboError_NoSpeedOrAcceleration
CWClassicFunctionNet, 677	CRoboDeviceNet, 485
ResetDacOffset	RoboError_OverPressure
COctoPotDeviceNet, 399	CRoboDeviceNet, 485
ResetHighpassFilter	RoboError ParameterNotAllowed
CFilterConfigurationNet, 138	CRoboDeviceNet, 485
ResetPipe	RoboError PeristalticTimeout
CGenericDevelopDeviceNet, 163	CRoboDeviceNet, 485
ResetStatus	RoboError Phase0OutOfRange
	CRoboDeviceNet, 485
CStg200xDownloadBasicNet, 563	Onodobeviceriet, 400

RoboError_PollLoopCanceled	SampleSize64Unsigned
CRoboDeviceNet, 486	Mcs::Usb, 81
RoboError_PollLoopCanceledAndStopMovement	SampleSizeNet
CRoboDeviceNet, 486	Mcs::Usb, 81
RoboError_Pressure	SBSVector1
CRoboDeviceNet, 486	Mcs::Usb, 58
RoboError_RangeExceeded	SBSVector2
CRoboDeviceNet, 486	Mcs::Usb, 58
RoboError_StateChangeNotPossible	SBSVector3
CRoboDeviceNet, 486	Mcs::Usb, 58
RoboError Timeout	SBSVector4
CRoboDeviceNet, 486	Mcs::Usb, 58
RoboError UnknownCommand	Scale
CRoboDeviceNet, 486	CCreateFilterNet, 118
Robolnject	ScanForHeadstages
Mcs::Usb, 79	CWClassicFunctionNet, 677
RoboMainLowLevelCommand	
	SCU1ElectrodeGroupHS1
CRoboDeviceNet, 487	Mcs::Usb, 82
RoboMainStatorLowLevelCommand	SCU1ElectrodeGroupHS2
CRoboStatorDeviceNet, 498	Mcs::Usb, 82
Robooycte2	SCU1ElectrodeGroupHS3
Mcs::Usb, 79	Mcs::Usb, 82
RoboStatusEvent	SCU1ElectrodeGroupHS4
CRoboDeviceNet, 487	Mcs::Usb, 82
RoboStatusEventDelegate	SCU1HS1Sideband1
Mcs::Usb, 95	Mcs::Usb, 83
Running	SCU1HS1Sideband1TriggerStatus
Mcs::Usb, 85	Mcs::Usb, 83
RunTable	SCU1HS1Sideband2
CRoboDacqNet, 465	Mcs::Usb, 83
1 /	SCU1HS1Sideband2TriggerStatus
s_FilterAttributesNet	Mcs::Usb. 83
s_FilterAttributesNet CFilterCoefficientsNet::s_FilterAttributesNet, 709	Mcs::Usb, 83 SCU1HS2Sideband1
-	SCU1HS2Sideband1
CFilterCoefficientsNet::s_FilterAttributesNet, 709	SCU1HS2Sideband1 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleDstSizeNet SampleRate	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1 SCU1HS3Sideband1TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1 SCU1HS3Sideband1TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2 Mcs::Usb, 83
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2 TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize32Signed Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2 Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize32Unsigned Mcs::Usb, 81 SampleSize32Unsigned Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2 Mcs::Usb, 83 SCU1HS4Sideband2 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2 Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus
CFilterCoefficientsNet::s_FilterAttributesNet, 709 SafeIS Mcs::Usb, 79 SampleDstSize16 Mcs::Usb, 81 SampleDstSize32 Mcs::Usb, 81 SampleDstSizeNet Mcs::Usb, 81 SampleRate CCreateFilterNet, 118 Samplerate CMcsUsbDacqNet, 297 SampleSize16Signed Mcs::Usb, 81 SampleSize16Unsigned Mcs::Usb, 81 SampleSize24Signed Mcs::Usb, 81 SampleSize24Unsigned Mcs::Usb, 81 SampleSize32Signed Mcs::Usb, 81 SampleSize32Signed Mcs::Usb, 81 SampleSize32Signed Mcs::Usb, 81	SCU1HS2Sideband1 Mcs::Usb, 83 SCU1HS2Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS2Sideband2 Mcs::Usb, 83 SCU1HS2Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband1 Mcs::Usb, 83 SCU1HS3Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS3Sideband2 Mcs::Usb, 83 SCU1HS3Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1 Mcs::Usb, 83 SCU1HS4Sideband1 TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband1TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2 Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus Mcs::Usb, 83 SCU1HS4Sideband2TriggerStatus

SCU1LEDStim1TriggerStatus	SCU2LEDStim2
Mcs::Usb, 83	Mcs::Usb, 84
SCU1LEDStim2	SCU2LEDStim2TriggerStatus
Mcs::Usb, 83	Mcs::Usb, 83
SCU1LEDStim2TriggerStatus	SCU2LEDStim3
Mcs::Usb, 83	Mcs::Usb, 84
SCU1LEDStim3	SCU2LEDStim3TriggerStatus
Mcs::Usb, 83	Mcs::Usb, 83
SCU1LEDStim3TriggerStatus	SCU2LEDStim4
Mcs::Usb, 83	Mcs::Usb, 84
SCU1LEDStim4	SCU2LEDStim4TriggerStatus
Mcs::Usb, 83	Mcs::Usb, 83
SCU1LEDStim4TriggerStatus	SCU_HeadstageIdEnumNet
Mcs::Usb, 83	Mcs::Usb, 81
SCU2ElectrodeGroupHS1	SCUDacqGroupChannelEnumNet
Mcs::Usb, 82	Mcs::Usb, 82
SCU2ElectrodeGroupHS2	SCUDigitalSourceEnumNet
Mcs::Usb, 82	Mcs::Usb, 82
SCU2ElectrodeGroupHS3	SegmentDefine
Mcs::Usb, 82	S .
•	CStg200xDownloadNet, 571
SCU2ElectrodeGroupHS4	SegmentSelect
Mcs::Usb, 82	CStg200xDownloadNet, 571
SCU2HS1Sideband1	SegmentStart
Mcs::Usb, 83	CStg200xDownloadNet, 572
SCU2HS1Sideband1TriggerStatus	SelectHeadstage
Mcs::Usb, 83	CW2100_FunctionNet, 620
SCU2HS1Sideband2	SelectTimeSlot
Mcs::Usb, 84	CW2100_StimulatorFunctionNet, 626
SCU2HS1Sideband2TriggerStatus	Send
Mcs::Usb, 83	CSerialPortNet, 520
SCU2HS2Sideband1	SendBuffered
Mcs::Usb, 84	CGilsonDeviceNet, 169
SCU2HS2Sideband1TriggerStatus	SendChannelData
Mcs::Usb, 83	CDigOutStimulatorFunctionNet, 131
SCU2HS2Sideband2	CStg200xDownloadBasicNet, 563
Mcs::Usb, 84	SendCommand
SCU2HS2Sideband2TriggerStatus	CLIH3DeviceNet, 201
Mcs::Usb, 83	SendImmediate
SCU2HS3Sideband1	CGilsonDeviceNet, 169
Mcs::Usb, 84	SendImmediateGetResponse
SCU2HS3Sideband1TriggerStatus	CGilsonDeviceNet, 169
Mcs::Usb, 83	SendMultiplexedData
SCU2HS3Sideband2	CStimulusFunctionNet, 581
Mcs::Usb, 84	SendPreparedData
SCU2HS3Sideband2TriggerStatus	CStimulusFunctionNet, 581
Mcs::Usb, 83	CW2100_StimulatorFunctionNet, 626
SCU2HS4Sideband1	SendStart
Mcs::Usb, 84	CStg200xBasicNet, 543
SCU2HS4Sideband1TriggerStatus	CStimulusFunctionNet, 581
Mcs::Usb, 83	CW2100_StimulatorFunctionNet, 627
SCU2HS4Sideband2	SendStartDacq
Mcs::Usb, 84	CMcsUsbDacqNet, 279
SCU2HS4Sideband2TriggerStatus	SendStartStgAndDacq
Mcs::Usb, 83	CMcsUsbDacqNet, 279
SCU2LEDStim1	SendStop
Mcs::Usb, 84	CStg200xBasicNet, 543
SCU2LEDStim1TriggerStatus	CStig200xBasicNet, 543 CStimulusFunctionNet, 582
Mcs::Usb. 83	CW2100 StimulatorFunctionNet. 627
17103030, 00	OVVETOU DITHUIATUT UTICITUTINEL DZ/

SendStopDacq	CGrapheneFunctionNet, 180
CMcsUsbDacqNet, 279	SetAdcOffset
SendStopStgAndDacq	CLIH3DeviceNet, 201
CMcsUsbDacqNet, 280	COctoPotDeviceNet, 399
SendStopStgAndDacqWithOptions	SetAdcOffsetPermanent
CMcsUsbDacqNet, 280	CLIH3DeviceNet, 202
SendSyncData	SetAdcSamplePos
CStg200xDownloadBasicNet, 564	CSafeISDeviceNet, 500
Sensor	SetAirpressureLimit
CFYIDeviceNet, 150	CRoboDeviceNet, 480
CMeasureTableDeviceNet, 370	SetAirValve
CPatchServerDeviceNet, 405	CRoboDeviceNet, 480
SerialNumber	SetAllDigout
CMcsUsbListEntryNet, 314	CRoboDacqNet, 465
CMcsUsbNet, 338	SetAmplificationSwitch
SerialPort	COctoPotDeviceNet, 399
CHLADeviceNet, 186	SetAmplitude
Set4ADCCatchampAverageShift	CChannelTestDeviceNet, 101
CMcsBus_SensorNet, 234	SetAmplitude_nA
Set4ADCMode	CTEERFunctionNet, 607
CMcsBus SensorNet, 234	SetAnalogOutADCRange
Set4DAC	CSCUFunctionNet, 516
CMcsBus_SensorNet, 234	SetAnalogOutChannel
Set_Values	CW2100_FunctionNet, 620
CNF_GenDeviceNet, 397	SetAnalogOutChannels
CPathIdentDeviceNet, 406	CSCUFunctionNet, 516
SetAbsMaxCurrentInMicroAmp	SetAnalogOutDACRange
CMultiwellOptoStimFunctionNet, 394	CSCUFunctionNet, 516
SetAccelerationI	SetAnalogOutFilter
CRoboStatorDeviceNet, 496	CW2100_FunctionNet, 620
SetAccelerationNativel	SetAnalogThresholdHigh
CRoboStatorDeviceNet, 496	CWarnerValveControllerDeviceNet, 662
SetAccelerationNativeXY	SetAnalogThresholdLow
CRoboStatorDeviceNet, 496	CWarnerValveControllerDeviceNet, 663
SetAccelerationNativeZ	SetAnalogVoltageRange
CRoboStatorDeviceNet, 496	CPPCFunctionNet, 431
SetAccelerationXY	, , , , , , , , , , , , , , , , , , ,
	SetAnalogVoltages
CRoboStatorDeviceNet, 496	CPPS_FunctionNet, 436
SetAccelerationZ	SetAttenuation
CRoboStatorDeviceNet, 496	CChannelTestDeviceNet, 101 SetAudioChannels
SetAccelGyroDesiredRate	
CW2100_FunctionNet, 620	CMeaAudioFunctionNet, 343
SetAccelGyroEnabled	CW2100_FunctionNet, 620
CW2100_FunctionNet, 620	SetAudioOutDacParameter
SetAccelRange	CLIH3DeviceNet, 202
CW2100_FunctionNet, 620	SetAutocalibrationDisabled
SetActiveRunningTableNumber	CStg200xBasicNet, 544
CWarnerValveControllerDeviceNet, 662	SetAxisConfig
SetADC	CRoboDeviceNet::RoboMainLowLevelCommands,
CWarnerValveControllerDeviceTesterFunctionNet,	706
673	SetAxisLED
SetAdcChannels	CRoboocyte2DeviceNet, 492
CSafeISDeviceNet, 500	SetAxisParametersEeprom
SetADCGain	CMcsBus_AxisParametersNet, 205, 206
CGrapheneFunctionNet, 179	SetBandwidthByIndex
SetADCInputOffset	CIntanMea_FunctionNet, 190
CCMOSMea_FunctionNet, 111	SetBaseFrequency
SetADCOffset	CRFFunctionNet, 452

SetBaseSamplerate	SetConfigurationBitBlu_Led
CCMOSMeaDeviceNet, 115	CRoboDacqNet, 466
CGrapheneASICDeviceNet, 171	SetConfigurationBitBlu_LedToggleFast
SetBath	CRoboDacqNet, 466
CCMOSMea_FunctionNet, 111	SetConfigurationBitBlu_LedToggleSlow
SetBathclamp	CRoboDacqNet, 466
COctoPotDeviceNet, 399	SetConfigurationBitCC_Gen
SetBathMode	CRoboDacqNet, 466
	SetConfigurationBitCV_Gen
CCMOSMea_FunctionNet, 111	
SetBlankingEnable	CRoboDacqNet, 466
CStg200xBasicNet, 544, 545	SetConfigurationBitRC_Gen
SetBuffer	CRoboDacqNet, 466
CGenericDevelopDeviceNet, 163	SetConfigurationBitRed_Led
SetBufferIndex	CRoboDacqNet, 466
CTEERFunctionNet, 607	SetConfigurationBitRed_LedSaturation
SetBusAddress	CRoboDacqNet, 466
CMcsBusNet, 245	SetConfigurationBitRed_LedToggleFast
SetBusAddressEeprom	CRoboDacqNet, 467
CMcsBusNet, 245	SetConfigurationBitRed_LedToggleSlow
SetByteBuffer	CRoboDacqNet, 467
CGenericDevelopDeviceNet, 163	SetConfigurationBitRelais
SetCalibration	CRoboDacqNet, 467
CTcxDeviceNet, 597	SetConfigurationBitRV_Gen
SetCardinalDacqSamplerate	CRoboDacqNet, 467
CInterfaceboardFunctionNet, 194	SetConfigurationBitStream
SetCardinalStgOutputrate	CRoboDacqNet, 467
CInterfaceboardFunctionNet, 194	SetConfigurationBitSupply
SetChannel	CRoboDacqNet, 467
CSw2to64DeviceNet, 585	SetControllerParams
	CTEERFunctionNet, 608
SetChannelmap	
CWClassicFunctionNet, 677	SetCrossTalkOffset
SetChannels	CRoboDacqNet, 467
CSw2to64DeviceNet, 586	SetCrossTalkOptimum
SetChannelSwitch	CRoboDacqNet, 467
COctoPotDeviceNet, 399	SetCurrentAirvalveLimit
SetChargingMode	CRoboDeviceNet, 480
CMultiBatteryChargerDeviceNet, 380	SetCurrentAndAir
SetChargingPCoefficient	CRoboDeviceNet, 480
CMultiBatteryChargerDeviceNet, 380	SetCurrentAndAirXY
SetCheckVoltage	CRoboStatorDeviceNet, 496
COkuvisionStimulatorDeviceNet, 403	SetCurrentEditTableNumber
SetClampMode	CWarnerValveControllerDeviceNet, 663
CTEERFunctionNet, 608	SetCurrentEnable
CWarnerUssingFunctionNet, 643	CTEERFunctionNet, 608
SetColorRgb	SetCurrentFactor
CMultiwellOptoStimFunctionNet, 395	COkuvisionStimulatorDeviceNet, 403
SetColorStr	SetCurrentMode
CMultiwellOptoStimFunctionNet, 395	CStg200xBasicNet, 545
•	
SetCommand	SetCurrentRangeSelectedIndex
CMcsBusNet, 246	CStg200xBasicNet, 546
CPedoterDeviceNet, 407	SetCycles
CRoboDacqNet, 465	CMeaCleanDeviceNet, 346
SetConfiguration	CMeaCoatDeviceNet, 351
CMcsUsbNet, 332	SetD
SetConfigurationBit	CTcxDeviceNet, 597
CRoboDacqNet, 465	SetDacAmplificationFactor
SetConfigurationBitAxc	CStg200xBasicNet, 546
CRoboDacqNet. 466	SetDacAutoControl

COsto Pat Daviso Not 200	CatDigitalDartDirection
COctoPotDeviceNet, 399 SetDacIdleValue	SetDigitalPortDirection
	CWarnerValveControllerDeviceNet, 664
CLIH3DeviceNet, 202	SetDigitalSource
SetDacMode	CMcsUsbDacqNet, 281–283
CSafeISDeviceNet, 500	SetDigitalStimulatorTrigger
SetDACOffset	CW2100_StimulatorFunctionNet, 627
CGrapheneFunctionNet, 180	SetDigitalStimulatorTriggerSlope
COkuvisionStimulatorDeviceNet, 403	CW2100_StimulatorFunctionNet, 627
SetDacOffset	SetDigout
CDacCalibrationFunctionNet, 120	CFluidControlDeviceNet, 147
CLIH3DeviceNet, 202	CRoboDacqNet, 467
COctoPotDeviceNet, 400	SetDigoutMode
SetDacOffsetPermanent	CStg200xBasicNet, 546
CLIH3DeviceNet, 203	SetDigOutState
SetDacPeriode	CLIH3DeviceNet, 203
CSafeISDeviceNet, 500	SetDigoutValue
SetDacPulseform	CStg200xBasicNet, 546
CSafeISDeviceNet, 500	SetDIO
SetDacqLegacyMode	CMcsBus_FYIExtensionNet, 208
CSCUFunctionNet, 516	SetDischargeCurrentSetPoint
SetDacRange	CMultiBatteryChargerDeviceNet, 380
CW2100_FunctionNet, 621	SetDisplayMode
SetDACs	CWarnerValveControllerDeviceNet, 664
CMcsBus_SensorNet, 234	SetDisplayText
SetDacUseIdleValue	CRoboDacqNet, 468
CLIH3DeviceNet, 203	SetDownsampleFactor
SetDacValue	CRoboDacqNet, 468
COctoPotDeviceNet, 400	SetDSPHighPassByIndex
SetDataMode	CIntanMea_FunctionNet, 191
CMcsUsbDacqNet, 280	SetDuration
SetDefault	CMeaCoatDeviceNet, 352
CWarnerValveControllerDeviceNet, 663	SetEEpromPage
SetDestinationSerialNumber	CLIH3DeviceNet, 203
CMcsUsbFactoryNet, 305	SetElectrodeDacMux
SetDetectionThreshold	CStg200xBasicNet, 547, 548
CMcsBus_SensorNet, 234	SetElectrodeEnable
SetDevice	CStg200xBasicNet, 549, 550
CTcxDeviceNet, 597	SetElectrodeMode
SetDeviceId	CStg200xBasicNet, 551, 552
CUsbDeviceConfigurationFunctionNet, 612	SetEnableAmplifierProtectionSwitch
SetDeviceList	CStg200xBasicNet, 553, 554
CPositionImpDeviceNet, 422	SetEnableHeaterLimit
SetDeviceName	CTcxDeviceNet, 597
CUsbDeviceConfigurationFunctionNet, 612	SetEnablePulse
SetDeviceType	CWarnerUssingFunctionNet, 643
CTcxDeviceNet, 597	SetEnableThermocouple
SetDevname	CTcxDeviceNet, 598
	SetExternalElectrodeEnable
CTcxDeviceNet, 597	
SetDiagnosticMode CIntanMea_FunctionNet, 190	CStg200xBasicNet, 554, 555 SetExternalLED
SetDigitalData CMooDigitalDataFunctionNet 262, 262	CTEERFunctionNet, 608
CMeaDigitalDataFunctionNet, 362, 363	SetFAAmplification
SetDigitalOut	CStg200xBasicNet, 555
CMeaDeviceNet, 357	SetFilter
SetDigitalOutPortInvert	CRoboDacqNet, 468
CWarnerValveControllerDeviceNet, 663	SetFilterCoeffs
SetDigitalOutPortValve	CRoboDacqNet, 468
CWarnerValveControllerDeviceNet, 664	SetFilterParameter

CFilterConfigurationNet, 138, 139	CMcsBusNet, 246
CFilterConfigurationRegisterNet, 140	SetHWSelectedChannels
SetFilterParameterPermanent	CWClassicFunctionNet, 678
CFilterConfigurationNet, 139	SetI
CFilterConfigurationRegisterNet, 140, 141	CTcxDeviceNet, 598
SetFilterParametersHeadstage	SetlClamp
CWClassicFunctionNet, 677	CRoboDacqNet, 468
SetFinalDischargeVoltage	SetICoeff
CMultiBatteryChargerDeviceNet, 381	CRobo FYITemp FunctionNet, 456
SetFrequency	SetICOffset
CChannelTestDeviceNet, 101	CRoboDacqNet, 468
CRadioControledDevicesNet, 445	SetIdleModeOffset
SetGain	CWarnerUssingFunctionNet, 644
CPgaDeviceNet, 411	SetIGain
SetGate	CRoboDacqNet, 468
CCMOSMea_FunctionNet, 111	SetImpedanceTestFrequency
SetGateFloating	CMealmpedanceDeviceNet, 369
CCMOSMea_FunctionNet, 111	SetImpld
SetGateToVOP	CPositionImpDeviceNet, 422
CCMOSMea_FunctionNet, 112	SetImplantCurrentSetpoint
SetGlobalRepeat	CPositionIIDeviceNet, 417
CDigOutStimulatorFunctionNet, 131	SetInMovement
SetGyroRange	CRoboDeviceNet, 480
CW2100_FunctionNet, 621	SetIntanRegister
SetHasChecksum	CIntanMea_FunctionNet, 191
CWClassicFunctionNet, 678	SetIntBuffer
SetHeadstage	CGenericDevelopDeviceNet, 164
CStg200xBasicNet, 555	SetIO
SetHeadstageFrameCyclesToComparePermanent	CWarnerValveControllerDeviceTesterFunctionNet,
CSCUFunctionNet, 517	673
SetHeadstageLinkSpeedPermanent	SetIODirection
CSCUFunctionNet, 517	CWarnerValveControllerDeviceTesterFunctionNet,
SetHeadstageNumberOfAnalogChannelsPermanent	674
CSCUFunctionNet, 517	SetloVoltage
SetHeadstageOnOff	CInterfaceboard2FunctionNet, 192
CW2100_FunctionNet, 621	SetLatency
CWClassicFunctionNet, 678	CMcsBus_SensorNet, 234
SetHeadstagePowerStateAtStart	SetLayoutConfiguration
CSCUFunctionNet, 517	CMEA2100x256FunctionNet, 340
SetHeadstageSampleratePermanent	SetLED
CSCUFunctionNet, 518	CRetinaLedDeviceNet, 447
SetHeadstageSamplingActive	SetLEDSwitch
CW2100_FunctionNet, 621	CMcsBus_ExtensionNet, 207
SetHeadstageToSleep	SetLength
CW2100_FunctionNet, 621	CRobo_FYIProgram_FunctionNet, 454
SetHeaterLimit	SetLiquidResistance
CTcxDeviceNet, 598	CTEERFunctionNet, 609
SetHighCurrentMode	CWarnerUssingFunctionNet, 644
CWarnerUssingFunctionNet, 643	SetListmodeIndexRange
SetHighpassFilterEnable	CStg200xBasicNet, 555
CFilterConfigurationNet, 139	SetListmodeTriggerSource
SetHWConfig	CStg200xBasicNet, 556
CRoboDeviceNet::RoboMainLowLevelCommands,	SetLowCurrentMode
706	CWarnerUssingFunctionNet, 644
SetHWRevision	SetLumi
CRoboDeviceNet::RoboMainLowLevelCommands,	CRetinaLedDeviceNet, 447
706	SetMaxCurrent
SetHWRevisionEeprom	CMeaCoatDeviceNet, 352

SetMaxDurationHighCurrentInMicroSec	CMcsBus_MotorControlNet, 221
CMultiwellOptoStimFunctionNet, 395	SetMCMaxCurrentEeprom
SetMaxDutyCycleHighCurrent	CMcsBus_MotorControlNet, 221
CMultiwellOptoStimFunctionNet, 395	SetMCMaxSpeed
SetMaxHeaterPowerMultiwell	CMcsBus_MotorControlNet, 222
CTcxDeviceNet, 598	SetMCMaxSpeedEeprom
SetMaxNoPressure	CMcsBus_MotorControlNet, 222
CRoboDeviceNet::RoboMainLowLevelCommands,	SetMCMaxTravel
706	CMcsBus_MotorControlNet, 222
SetMaxNoPressureWaitTime	SetMCMaxTravelEeprom
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus MotorControlNet, 222
706	SetMCMaxTravelShortCommand
SetMaxP	CMcsBus_MotorControlNet, 222
CTcxDeviceNet, 598	SetMCNewPosition
SetMaxPower	CMcsBus_MotorControlNet, 222
COkuvisionStimulatorDeviceNet, 404	SetMCOutputOnOff
CRobo_FYITemp_FunctionNet, 456	CMcsBus_MotorControlNet, 223
SetMaxPressureWaitTime	SetMCReference
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus MotorControlNet, 223
706	SetMCReferenceCurrent
SetMaxVoltage	CMcsBus_MotorControlNet, 223
CMeaCleanDeviceNet, 347	SetMCReferenceCurrentEeprom
COkuvisionStimulatorDeviceNet, 404	CMcsBus_MotorControlNet, 223
SetMCAcceleration	SetMCRegulatorGain
CMcsBus_MotorControlNet, 218	CMcsBus_MotorControlNet, 223
SetMCAccelerationEeprom	SetMCRegulatorGainEeprom
CMcsBus_MotorControlNet, 219	CMcsBus_MotorControlNet, 223
SetMCAccelerationShortCommand	SetMCRotation
CMcsBus_MotorControlNet, 219	CMcsBus_MotorControlNet, 224
SetMCAxisRevisionEeprom	SetMCScalingFactor
CMcsBus_MotorControlNet, 219	CMcsBus_MotorControlNet, 224
SetMCBreakCurrent	SetMCScalingFactorEeprom
CMcsBus_MotorControlNet, 219	CMcsBus MotorControlNet, 224
SetMCBreakCurrentEeprom	SetMCSpeed
CMcsBus_MotorControlNet, 219	CMcsBus_MotorControlNet, 224
SetMCConfig	SetMCSpeedEeprom
_	CMcsBus_MotorControlNet, 224
CMcsBus_MotorControlNet, 219	
SetMCConfigEeprom CMcsBus MotorControlNet, 220	SetMCSpeedShortCommand CMcsBus MotorControlNet, 224
	——————————————————————————————————————
SetMCCurrent CMooRya MotorControlNet 220	SetMCSpeedUnitEeprom
CMcsBus_MotorControlNet, 220	CMcsBus_MotorControlNet, 225
SetMCCurrentEeprom CMosRys MaterCentralNet 220	SetMCStandbyCurrent CMacRus MeterControlNet 225
CMcsBus_MotorControlNet, 220	CMcsBus_MotorControlNet, 225
SetMCCurrentMode CMacRus MaterControlNet 000	SetMCStandbyCurrentEeprom
CMcsBus_MotorControlNet, 220	CMcsBus_MotorControlNet, 225
SetMCCurrentModeEeprom	SetMCStandbyTime
CMcsBus_MotorControlNet, 220	CMcsBus_MotorControlNet, 225
SetMCCurrentModeShortCommand	SetMCStandbyTimeEeprom
CMcsBus_MotorControlNet, 220	CMcsBus_MotorControlNet, 225
SetMCCurrentPosition	SetMeasurementMode
CMcsBus_MotorControlNet, 221	CStg200xBasicNet, 557
SetMCCurrentShortCommand	SetMinimalThreshold
CMcsBus_MotorControlNet, 221	CMcsBus_SensorNet, 234
SetMCMaxAcceleration	SetMinNoPressureWaitTime
CMcsBus_MotorControlNet, 221	CRoboDeviceNet::RoboMainLowLevelCommands,
SetMCMaxAccelerationEeprom	706
CMcsBus_MotorControlNet, 221	SetMinPressure
SetMCMaxCurrent	CRoboDeviceNet, 481

CRoboDeviceNet::RoboMainLowLevelCommands,	SetPidParameter
706	COctoPotDeviceNet, 400
SetMinPressureWaitTime	SetPiezoState
CRoboDeviceNet::RoboMainLowLevelCommands,	CMcsBus_SensorNet, 235
706	SetPlateMux
SetMinVoltage	CMultiwellDeviceNet, 389
CMeaCleanDeviceNet, 347	SetPlateType
SetModeSelect	CMultiwellDeviceNet, 390
CPulseGeneratorFunctionNet, 443	SetPoti
SetMovePump	CMcsUsbDacqNet, 283
CMcsBus SensorNet, 235	SetPowerMuxPlate
SetMultiHeadstageMode	CMultiwellDeviceNet, 390
CW2100 FunctionNet, 621	SetPowerStrength
SetNanoVoltsPerKelvin	CPositionIIDeviceNet, 417
CMcsBus_TempSensorNet, 238	SetPressureOffset
SetNeurochipMemoryData	CMcsBus_SensorNet, 235
CCMOSMea_FunctionNet, 112	CPPCFunctionNet, 431
SetNoFilterCoeffs	SetPressureRange
CRoboDacqNet, 468	CPPCFunctionNet, 431
SetNumberOfAnalogChannels	SetPulse
CMeaDeviceNet, 358	CWarnerUssingFunctionNet, 644
SetNumberOfChannels	SetPulseform
CMeaDeviceNet, 359	COkuvisionStimulatorDeviceNet, 404
COctoPotDeviceNet, 400	SetPulseLength
SetOffsetCurrent	CPulseGeneratorFunctionNet, 443
CMeaCoatDeviceNet, 352	SetPumpCouple
SetOnOff	CPPS_FunctionNet, 437
CTcxDeviceNet, 598	SetPumpEnableSpeedRatio
SetOutputMap	CPPS_FunctionNet, 437
CStg200xDownloadNet, 572	SetPumpFastOnOff
SetOutputRate	CPPS_FunctionNet, 437
COctoPotDeviceNet, 400	SetPumpFastSpeed
CStg200xBasicNet, 557	CPPS_FunctionNet, 437
SetP	SetPumpFunctionSpeeds
CTcxDeviceNet, 599	CPPS_FunctionNet, 437
SetParameter	SetPumpManualOnOff
CRoboDeviceNet::RoboMainLowLevelCommands,	CPPS_FunctionNet, 437
707	SetPumpMaxSpeed
SetPattern	CPPS_FunctionNet, 437
	SetPumpModeType
CMeaSwitchDeviceNet, 371 SetPatternBool	CPPCFunctionNet, 432
CMeaSwitchDeviceNet, 372	CPPS FunctionNet, 437
SetPatternListEntry	SetPumpSpeed
COctoPotDeviceNet, 400	CRoboFluidDeviceNet, 489
SetPauseDuration	SetPumpSpeedRatio
CMeaCoatDeviceNet, 352	CPPS_FunctionNet, 438
SetPCoeff	SetPumpSpeedUnit
CRobo_FYITemp_FunctionNet, 456	CPPCFunctionNet, 432
SetPeriod	CPPS_FunctionNet, 438
CPulseGeneratorFunctionNet, 443	SetPWM
SetPeriod_us	CFluidControlDeviceNet, 147
CTEERFunctionNet, 609	SetRampParameter
SetPermanentCurrentInMicroAmp	COctoPotDeviceNet, 400
CMultiwellOptoStimFunctionNet, 395	SetRatedCapacity
SetPersistency	CMultiBatteryChargerDeviceNet, 381
CRetinaLedDeviceNet, 447	SetRatedCapacityVolatile
SetPGain	CMultiBatteryChargerDeviceNet, 381
CRoboDacqNet, 468	SetRecordingNumber

CRoboDacqNet, 469	SetSearchReferenceMethod
SetReferenceElectrodeMode	CRoboDeviceNet::RoboMainLowLevelCommands
CSCUFunctionNet, 518	707
SetReferenceElectrodeSwitchState	SetSearchReferenceMoveOut
CSCUFunctionNet, 518	CRoboDeviceNet::RoboMainLowLevelCommands
SetRegionOfInterests	708
CCMOSMeaDeviceNet, 115	SetSearchReferenceOffsetPos
	CRoboDeviceNet::RoboMainLowLevelCommands
CGrapheneASICDeviceNet, 171	
SetRegulationTimeouts	708
CMcsBus_SensorNet, 235	SetSelectedChannels
SetRegulatorFactor	CMcsUsbDacqNet, 285–287
CMcsBus_SensorNet, 235	CW2100_FunctionNet, 621
SetRegulatorOnOff	SetSelectedChannelsQueue
CMcsBus_SensorNet, 235	CMcsUsbDacqNet, 287–289
CRobo_FYITemp_FunctionNet, 456	SetSelectedData
SetRepeat	CMcsUsbDacqNet, 289-291
CRetinaLedDeviceNet, 448	SetSelectedHeadstage
SetRepeats	CWClassicFunctionNet, 679
CProgramPressureCurveNet, 440	SetSensorType
SetResetFilter	CTcxDeviceNet, 599
CWClassicFunctionNet, 678	SetSerialNumberHeadstage
SetRFFrequency	CWClassicFunctionNet, 679
CPositionImpDeviceNet, 423	SetSetpoint
SetRFFrequencyHeadstage	CTcxDeviceNet, 599
CWClassicFunctionNet, 678	SetShortBuffer
SetRFFrequencyReceiver	CGenericDevelopDeviceNet, 164
CWClassicFunctionNet, 678	SetSimulation
SetRFFrequencyReceiverEeprom	CRoboDacqNet, 469
CWClassicFunctionNet, 678	SetSineParameter
SetRFLostBehaviour	COctoPotDeviceNet, 400
CWClassicFunctionNet, 678	SetSingleHeater
SetRFPower	CMcsBus_FYIExtensionNet, 208
CWClassicFunctionNet, 679	SetSingleValve
SetRotatePump	CFluidControlDeviceNet, 148
CMcsBus_SensorNet, 235	CRoboFluidDeviceNet, 489
SetRTC	SetSlope
COkuvisionStimulatorDeviceNet, 404	CMeaCleanDeviceNet, 347
CPositionIIDeviceNet, 417	CMeaCoatDeviceNet, 353
SetSampleInterval	SetSoftwareKey
CLIH3DeviceNet, 204	CMcsUsbNet, 332
SetSamplePeriode	SetSollPressure
CMcsBus_SensorNet, 236	CMcsBus_SensorNet, 236
SetSamplerate	SetSollTemp
CMcsUsbDacqNet, 283	CRobo_FYITemp_FunctionNet, 456
SetScreen	SetSourceBulk
	CCMOSMea_FunctionNet, 112
CRoboDacqNet, 469	
SetSearchReferenceFastAccel	SetSourceDrain
CRoboDeviceNet::RoboMainLowLevelCommands,	CCMOSMea_FunctionNet, 112
707	SetSourceGate
SetSearchReferenceFastSpeed	CCMOSMea_FunctionNet, 112
CRoboDeviceNet::RoboMainLowLevelCommands,	SetSpeedI
707	CRoboStatorDeviceNet, 497
SetSearchReferenceFineAccel	SetSpeedNativel
CRoboDeviceNet::RoboMainLowLevelCommands,	CRoboStatorDeviceNet, 497
707	SetSpeedNativeXY
SetSearchReferenceFineSpeed	CRoboStatorDeviceNet, 497
CRoboDeviceNet::RoboMainLowLevelCommands,	SetSpeedNativeZ
707	CRoboStatorDeviceNet. 497
101	OI WOOGAWA DEVICEINEL 43/

SetSpeedXY	CWarnerValveControllerDeviceTesterFunctionNet,
CRoboStatorDeviceNet, 497	674
SetSpeedZ	SetUByteBuffer
CRoboStatorDeviceNet, 497	CGenericDevelopDeviceNet, 165
SetStartTriggerSlope	SetUClamp
CDigOutStimulatorFunctionNet, 131	CRoboDacqNet, 469
SetStateDebugData	SetUCOffset
CPositionIIDeviceNet, 419	CRoboDacqNet, 469
SetStateEventData	SetUIntBuffer
CPositionIIDeviceNet, 419	CGenericDevelopDeviceNet, 165
SetStgProgramInfo	SetupGroupDacqQueue
CStg200xBasicNet, 557	CMcsUsbDacqNet, 291
SetStimulusSites	SetupRetriggerMode
CCMOSMea_FunctionNet, 112	CStg200xDownloadBasicNet, 564
SetStopTriggerSlope	SetupTrigger
CDigOutStimulatorFunctionNet, 132	CStg200xDownloadBasicNet, 565
SetStringFormat	CStimulusFunctionNet, 582
CMcsUsbListEntryNet, 313	SetupTriggerSingle
CMcsUsbListNet, 317	CStg200xDownloadBasicNet, 565
SetSubChannel	CStimulusFunctionNet, 583
CMcsBus_MotorControlNet, 225	SetUseBubble
SetSwitches	CPPS FunctionNet, 438
CSafeISDeviceNet, 501	SetUserParameter
SetSyncoutMap	CRoboDeviceNet::RoboMainLowLevelCommands,
CStg200xBasicNet, 557	708
SetTableName	SetUShortBuffer
CWarnerValveControllerDeviceNet, 664	CGenericDevelopDeviceNet, 166
SetTablepointer	SetUVOffset
CRetinaLedDeviceNet, 448	CRoboDacqNet, 469
SetTableStep	SetValue
CWarnerValveControllerDeviceNet, 665	
	CGenericDevelopDeviceNet, 167
SetTableStepAll	SetValve
CWarnerValveControllerDeviceNet, 665	CFluidControlDeviceNet, 148
SetTestMode	CRoboFluidDeviceNet, 489
CRFFunctionNet, 452	SetValve1
SetThermocoupleNanovoltPerKelvin	CRobo_FYIProgram_FunctionNet, 454
CFluidControlDeviceNet, 148	SetValve2
CTcxDeviceNet, 599	CRobo_FYIProgram_FunctionNet, 454
SetThermoOffset	SetValveActive
CMcsBus_TempSensorNet, 238	CPPCFunctionNet, 432
Settings	CWarnerValveControllerDeviceNet, 665
Mcs::Usb, 94	SetValveCurrent
SetTouchPadEnable	CWarnerValveControllerDeviceNet, 665
CMultiwellDeviceNet, 390	SetValveDigitalInInvert
SetTransformer	CWarnerValveControllerDeviceNet, 666
CMeFunctionNet, 374	SetValveDigitalInPort
SetTrigger	CWarnerValveControllerDeviceNet, 666
CRetinaLedDeviceNet, 448	SetValveLedOn
CWarnerValveControllerDeviceTesterFunctionNet,	CWarnerValveControllerDeviceNet, 666
674	SetValveManualGroup
SetTriggerMaskValue	CWarnerValveControllerDeviceNet, 666
CMeaDeviceNet, 359	SetValveManualState
CRoboDacqNet, 469	CWarnerValveControllerDeviceNet, 667
SetTriggerPeriod	SetValveMode
CMeaDeviceNet, 360	CWarnerValveControllerDeviceNet, 667
SetTriggerSource	SetValves
CStg200xBasicNet, 558	CMcsBus_FYIExtensionNet, 208
SetTriggerSyncDirection	SetValvesActiveMap
	COLTAITOUT IOLITOITIAD

CWarnerValveControllerDeviceNet, 667	CMeaUSBDeviceNet, 373
SetValvesManualStateMap	SetVoltageRangeInMicroVolt
CWarnerValveControllerDeviceNet, 667	CMcsUsbDacqNet, 292
SetValveTableEntry	SetVoltageRangeSelectedIndex
CWarnerValveControllerDeviceNet, 667	CStg200xBasicNet, 558
SetVdsVgs	SetVoltageResolution
CGrapheneFunctionNet, 181	CGrapheneFunctionNet, 183
SetVdVs	SetVoltageRs485ALimit
CGrapheneFunctionNet, 181	CRoboDeviceNet, 481
SetVdVsDAC	SetVoltageRs485BLimit
CGrapheneFunctionNet, 182	CRoboDeviceNet, 481
SetVelocityI	SetVoltageValvesLimit
CRoboStatorDeviceNet, 497	CRoboDeviceNet, 481
SetVelocityXY	SetWaveform
CRoboStatorDeviceNet, 497	CChannelTestDeviceNet, 101
SetVelocityZ	CTEERFunctionNet, 609
CRoboStatorDeviceNet, 497	SetWaveLengthInNanometer
SetVMMaxNegativeCurrent	CMultiwellOptoStimFunctionNet, 396
CMcsBus_VoltageModeNet, 241	SetWorkingFrequency
SetVMMaxNegativeCurrentEeprom	CRFFunctionNet, 452
CMcsBus_VoltageModeNet, 241	SetWPADebugMode
SetVMMaxNegativeVoltage	CWClassicFunctionNet, 679
CMcsBus_VoltageModeNet, 241	SetWPAType
SetVMMaxNegativeVoltageEeprom	CWClassicFunctionNet, 679
CMcsBus_VoltageModeNet, 241	SetXGain
SetVMMaxPositiveCurrent	CRoboDacqNet, 469
CMcsBus_VoltageModeNet, 242	Sideband
SetVMMaxPositiveCurrentEeprom	CStimulusFunctionNet::SidebandData, 711
CMcsBus_VoltageModeNet, 242	SidebandData
SetVMMaxPositiveVoltage	CStimulusFunctionNet::SidebandData, 711
CMcsBus_VoltageModeNet, 242	Signed_16bit
SetVMMaxPositiveVoltageEeprom	Mcs::Usb, 61
CMcsBus VoltageModeNet, 242	Signed 24bit
SetVMOutputOnOff	Mcs::Usb, 61
CMcsBus_VoltageModeNet, 242	Signed_32bit
SetVMVoltage	Mcs::Usb, 61
CMcsBus_VoltageModeNet, 242	Sine
SetVolatileClampOffset	Mcs::Usb, 88
CMultiwellDeviceNet, 391	SineStart
SetVoltage12VLimit	COctoPotDeviceNet, 401
CRoboDeviceNet, 481	Single
SetVoltage5VLimit	Mcs::Usb, 80
CRoboDeviceNet, 481	SingleWell
SetVoltageAirvalveLimit	Mcs::Usb, 86
CRoboDeviceNet, 481	SixWell
SetVoltageClampControllerParam D	
_	Mcs::Usb, 86
CWarnerUssingFunctionNet, 645	SiZe
SetVoltageClampControllerParam_I	DigitalSource< digitalsourceenum >, 682
CWarnerUssingFunctionNet, 645	DigitalSourceGeneral, 684
SetVoltageClampControllerParam_P	SmartImplant
CWarnerUssingFunctionNet, 645	Mcs::Usb, 79 SN
SetVoltageMode	_
CStg200xBasicNet, 558	HeadStageIDType, 696
SetVoltageRange	SOFAndCTRLword
CGrapheneFunctionNet, 182, 183	Mcs::Usb, 60
SetVoltageRangeByIndex	Software Macrillah 66
CMcsUsbDacqNet, 291	Mcs::Usb, 66
SetVoltageRangeForIndex	SoftwareDongle

Maguellah 70	CMaal labNat 205
Mcs::Usb, 79	CMcsUsbNet, 335
Source	Status_DataUnderrun
DigitalSource< digitalsourceenum >, 683	CMcsUsbNet, 335
DigitalSourceGeneral, 684	Status_DeviceLocked
Standby	CMcsUsbNet, 335
Mcs::Usb, 81, 90	Status_DeviceNotFound
Start CMacClass Pavias Nat. 847	CMcsUsbNet, 335
CMeaCleanDeviceNet, 347	Status_DeviceRemoved
CMeaCoatDeviceNet, 353	CMcsUsbNet, 335
CRobo_FYIProgram_FunctionNet, 454	Status_DevNotResponding
Mcs::Usb, 63	CMcsUsbNet, 335
StartDacq CMcsUsbDacqNet, 292, 293	Status_EndpointHalted CMcsUsbNet, 336
StartInternalCalibration	Status_ErrorBusy
CTEERFunctionNet, 609	CMcsUsbNet, 336
StartLoop	Status ErrorShortTransfer
CMcsUsbDacqNet, 294, 295	CMcsUsbNet, 336
StartMCMovement	Status Fifo
CMcsBus_MotorControlNet, 226	CMcsUsbNet, 336
StartMeasurement	Status FrameControlOwned
CMealmpedanceDeviceNet, 369	CMcsUsbNet, 336
StartPoll	Status InternalHcError
CStimulusFunctionNet, 583	CMcsUsbNet, 336
CW2100 StimulatorFunctionNet, 627	Status InvalidDeviceHandle
StartQueue	CMcsUsbNet, 336
CRoboDeviceNet, 481	Status InvalidHandle
StartSampling	CMcsUsbNet, 336
CTEERFunctionNet, 609	Status InvalidParameter
StartSync	CMcsUsbNet, 336
CMcsBus_SensorNet, 236	Status_InvalidPipeHandle
State	CMcsUsbNet, 336
HeadStageIDTypeState, 698	Status_InvalidUrbFunction
Mcs::Usb, 91	CMcsUsbNet, 337
Status	Status_IoPending
CUsbExceptionNet, 613	CMcsUsbNet, 337
Status_AlreadyConfigured	Status_IoTimeout
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_BadStartFrame	Status_IsochRequestFailed
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_Btstuff	Status_LastUsbErrorMismatch
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_BufferOverrun	Status_NoBandwidth
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_BufferUnderrun	Status_NoMemory
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_Canceled	Status_NoSuchDevice
CMcsUsbNet, 334	CMcsUsbNet, 337
Status_Canceling	Status_NotAccessed
CMcsUsbNet, 335	CMcsUsbNet, 337
Status_ConnectedPipes	Status_NotSupported
CMcsUsbNet, 335	CMcsUsbNet, 337
Status_ControlNotOwned	Status_PidCheckFailure
CMcsUsbNet, 335	CMcsUsbNet, 337
Status_Crc	Status_PipeNotLinked
CMcsUsbNet, 335	CMcsUsbNet, 338
Status_DataOverrun	Status_RequestFailed
CMcsUsbNet, 335	CMcsUsbNet, 338
Status_DataToggleMismatch	Status_RequestMutexFailed

CMcsUsbNet, 338	Mcs::Usb, 77
Status_RequestMutexTimeout	STG5
CMcsUsbNet, 338	Mcs::Usb, 77
Status_Stall	STG_DestinationEnumNet
CMcsUsbNet, 338	Mcs::Usb, 85
Status Unconfigured	StgListModeTrigger
CMcsUsbNet, 338	Mcs::Usb, 64
Status_UnexpectedPid	StgStatusNet, 712
CMcsUsbNet, 338	FromIntPtr, 712
STG	FromPtr, 712
Mcs::Usb, 77	ListOfChangedTriggers, 712
Stg1	TiggerStatus, 712
Mcs::Usb, 64	StgTrigger
STG1DACSignalGroup	Mcs::Usb, 64
Mcs::Usb, 61, 71, 82	StillConnected
STG1SidebandsGroup	CRadioControledDevicesNet, 445
Mcs::Usb, 61, 71, 82	Stimulation
STG1TriggerStatusGroup	Mcs::Usb, 57
Mcs::Usb, 61, 71, 82	StimulationLayoutConfigurationEnumNet
Stg2	Mcs::Usb, 86
Mcs::Usb, 64	Stimulator
Stg200xDigoutModeEnumNet	CW2100_FunctionNet, 622
Mcs::Usb, 84	Stimulus
Stg200xMultiFileSubmodeEnumNet	CCMOSMeaDeviceNet, 116
Mcs::Usb, 84	CStg200xDownloadBasicNet, 566
Stg200xPollStatusEvent	Stimulus Device Data And Unrolled Data
CStg200xDownloadNet, 573	CStimulusFunctionNet::StimulusDeviceDataAndUnrolledData,
Stg200xSegmentFlagsEnumNet	713 StimulusFunction
Mcs::Usb, 84	
Stg200xTriggerStatusEnumNet Mcs::Usb, 85	CLIH3DeviceNet, 204 StimulusParameters
STG2DACSignalGroup	HeadStageIDType, 696
Mcs::Usb, 71, 82	Stop
STG2SidebandsGroup	CMeaCleanDeviceNet, 347
Mcs::Usb, 71, 82	CMeaCoatDeviceNet, 353
STG2TriggerStatusGroup	Mcs::Usb, 63, 76, 80
Mcs::Usb, 71, 82	StopDacq
Stg3	CMcsUsbDacqNet, 296
Mcs::Usb, 64	StopLoop
STG3008 FA	CMcsUsbDacqNet, 296
Mcs::Usb, 77	StopMCMovement
Stg3008FilterAmpAmplificationEnumNet	CMcsBus_MotorControlNet, 226
Mcs::Usb, 85	StopMovement
STG4002	CRoboDeviceNet, 482
Mcs::Usb, 77	StopMovementI
STG4002_opto	CRoboStatorDeviceNet, 498
Mcs::Usb, 77	StopMovementXY
STG4004	CRoboStatorDeviceNet, 498
Mcs::Usb, 77	StopMovementZ
STG4004_opto	CRoboStatorDeviceNet, 498
Mcs::Usb, 77	StopPlateClamp
STG4008	CMultiwellDeviceNet, 391
Mcs::Usb, 77	StopPoll
STG4008_opto	CStimulusFunctionNet, 583
Mcs::Usb, 77	CW2100_StimulatorFunctionNet, 627
STG400x	StopSampling
Mcs::Usb, 77	CTEERFunctionNet, 610
STG400x_opto	StopTable

CRoboDacqNet, 469, 470	Mcs::Usb, 52
StorageCharge	TBSI_15
Mcs::Usb, 70	Mcs::Usb, 52
StoreValveTable	TBSI_31
CWarnerValveControllerDeviceNet, 668	Mcs::Usb, 52
SubtractFromAll	TBSI_5
Mcs::Usb, 80	Mcs::Usb, 52
SubtractFromAllOther	TBSI_63
Mcs::Usb, 80	Mcs::Usb, 52
SubtractFromReferenceElectrodeOnly	TBSI_Dacq
Mcs::Usb, 80	Mcs::Usb, 78
SubtractionOff	TBSI_DACQDigitalSourceEnumNet
Mcs::Usb, 80	Mcs::Usb, 86
SuperSpeed	TBSI_Reserved
Mcs::Usb, 70	Mcs::Usb, 52
Sw2to64	TbsiDacq
Mcs::Usb, 78	Mcs::Usb, 68
SwitchOnOff	TbsiDacqHeadstage
CPositionIIDeviceNet, 419	Mcs::Usb, 67
SYNC_BIT0	TbsiDacqInterfaceboard
CW2100_StimulatorFunctionNet, 628	Mcs::Usb, 67
SYNC_BIT1	TC01
CW2100_StimulatorFunctionNet, 628	Mcs::Usb, 77
SYNCOUT1	TC02
Mcs::Usb, 84	Mcs::Usb, 77
SYNCOUT2	TCX
Mcs::Usb, 84	Mcs::Usb, 77
SYNCOUT3	TcxDeviceTypeEnumNet
Mcs::Usb, 84	Mcs::Usb, 87
SYNCOUT4	TcxSensorTypeEnumNet
Mcs::Usb, 84	Mcs::Usb, 87
SYNCOUT5	TeerClampModeEnumNet
Mcs::Usb, 84	Mcs::Usb, 87
SYNCOUT6	TEERFunctionNet
Mcs::Usb, 84	CTEERMachineDeviceNet, 611
SYNCOUT7	TeerWaveformEnumNet
Mcs::Usb, 84	Mcs::Usb, 87
SYNCOUT8	Tersens
Mcs::Usb, 84	Mcs::Usb, 77
syncoutdata	Test_ADC_EPC10
Mcs::Usb, 85	Mcs::Usb, 68
SyncStart	Test_DAC_EPC10
Mcs::Usb, 85	Mcs::Usb, 69
	ThrowCUsbExceptionNetOnError
Table	CMcsUsbFunctionNet, 309
Mcs::Usb, 94	CMcsUsbNet, 332
Table_Wait	TiggerStatus
CRoboDacqNet, 470	StgStatusNet, 712
TableDefBegin	TimeResolutionInNanoSeconds
CRoboDacqNet, 470	W2100_StimulusParametersNet, 715
TableDefEnd	Timestamp
CRoboDacqNet, 470	Mcs::Usb, 60
TableEntryChangedEvent	ТоСрр
CWarnerValveControllerDeviceNet, 671	CFilterCoefficientsNet::s_FilterAttributesNet, 709
TactSwitchGetState	ToString
CMcsBus_SensorNet, 236	CFilterPropertyNet, 142
TactSwitchSetDisplay	CMcsUsbListEntryNet, 314
CMcsBus_SensorNet, 236	HeadStageIDType, 695
TBSI_127	O 715-7

HeadstageIDTypeObject, 697	tsDigitalIn16
TouchTest	Mcs::Usb, 88
Mcs::Usb, 94	tsDigitalIn17
Triggerbox_AMS	Mcs::Usb, 88
Mcs::Usb, 78	tsDigitalIn18
Triggerbox_AMS3	Mcs::Usb, 88
Mcs::Usb, 78	tsDigitalIn19
Triggerbox_IMS	Mcs::Usb, 88
Mcs::Usb, 78	tsDigitalIn2
Triggerbox_R5	Mcs::Usb, 88
Mcs::Usb, 78	tsDigitalIn20
TriggerMask Default	Mcs::Usb, 88
CRoboDacqNet, 470	tsDigitalIn21
TriggerOnly	Mcs::Usb, 88
Mcs::Usb, 85	tsDigitalIn22
TriggerSourceEnumNet	Mcs::Usb, 88
Mcs::Usb, 88	tsDigitalIn23
TriggerStatus	Mcs::Usb, 88
CMcsUsbDeviceStatePushFunctionNet, 298	tsDigitalIn24
CMcsUsbDeviceStatePushNet, 299	Mcs::Usb, 88
TriggerStatus1	tsDigitalIn25
Mcs::Usb, 59	Mcs::Usb, 88
TriggerStatus2	tsDigitalIn26
Mcs::Usb, 59	Mcs::Usb, 88
TriggerStatus3	tsDigitalIn27
Mcs::Usb, 59	Mcs::Usb, 88
TriggerStatus4	tsDigitalIn28
Mcs::Usb, 59	Mcs::Usb, 88
TriggerValue_MoveAbs	tsDigitalIn29
CRoboDacqNet, 470	Mcs::Usb, 88
TriggerValue_StartQueue	tsDigitalIn3
CRoboDacqNet, 470	Mcs::Usb, 88
tsAuxIn1	tsDigitalIn30
Mcs::Usb, 89	Mcs::Usb, 88
tsAuxIn2	tsDigitalIn31
Mcs::Usb, 89	Mcs::Usb, 88
tsDACQCy1Dev1Runs	tsDigitalIn32
Mcs::Usb, 90	Mcs::Usb, 88
tsDACQCy1Dev2Runs	tsDigitalIn4
Mcs::Usb, 90	Mcs::Usb, 88
tsDACQCy2Dev1Runs	tsDigitalIn5
Mcs::Usb, 90	Mcs::Usb, 88
tsDACQCy2Dev2Runs	tsDigitalIn6
Mcs::Usb, 90	Mcs::Usb, 88
tsDigitalIn1	tsDigitalIn7
Mcs::Usb, 88	Mcs::Usb, 88
tsDigitalIn10	tsDigitalIn8
Mcs::Usb, 88	Mcs::Usb, 88
tsDigitalIn11	tsDigitalIn9
Mcs::Usb, 88	Mcs::Usb, 88
tsDigitalIn12	tsDigitalPuse0
Mcs::Usb, 88	Mcs::Usb, 89
tsDigitalIn13	tsDigitalPuse1
Mcs::Usb, 88	Mcs::Usb, 89
tsDigitalIn14	tsDigitalPuse10
Mcs::Usb, 88	Mcs::Usb, 89
tsDigitalIn15	tsDigitalPuse11
Mcs::Usb, 88	Mcs::Usb, 89

1 D: 1: ID	. = " 140
tsDigitalPuse12	tsFeedback10
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse13	tsFeedback11
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse14	tsFeedback12
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse15	tsFeedback13
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse16	tsFeedback14
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse17	tsFeedback15
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse18	tsFeedback16
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse19	tsFeedback17
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse2	tsFeedback18
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse20	tsFeedback19
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse21	tsFeedback2
Mcs::Usb, 90	Mcs::Usb, 88
tsDigitalPuse22	tsFeedback20
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse23	tsFeedback21
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse24	tsFeedback22
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse25	tsFeedback23
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse26	tsFeedback24
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse27	tsFeedback25
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse28	tsFeedback26
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse29	tsFeedback27
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse3	tsFeedback28
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse30	tsFeedback29
Mcs::Usb, 90	Mcs::Usb, 89
tsDigitalPuse31	tsFeedback3
Mcs::Usb, 90	Mcs::Usb, 88
tsDigitalPuse4	tsFeedback30
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse5	tsFeedback31
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse6	tsFeedback32
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse7	tsFeedback4
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse8	tsFeedback5
Mcs::Usb, 89	Mcs::Usb, 89
tsDigitalPuse9	tsFeedback6
Mcs::Usb, 89	Mcs::Usb, 89
tsFeedback1	tsFeedback7
Mcs::Usb, 88	Mcs::Usb, 89
IVIUSUSD, 00	IVICSUSD, 89

tsFeedback8	UpdateTrigger
Mcs::Usb, 89	Mcs::Usb, 85
tsFeedback9	USB
Mcs::Usb, 89	FirmwareDestinationNames, 693
tsNone	Mcs::Usb, 53
Mcs::Usb, 88	USB_TARGET1
tsSidebandBit8	Mcs::Usb, 57
Mcs::Usb, 90	USB_TARGET2
tsTriggered	Mcs::Usb, 57
Mcs::Usb, 90	USB_TARGET3
TxnGetSerialNumber	Mcs::Usb, 57
CMcsUsbNet, 332	usbSetupPacket_t, 714
TxnSetSerialNumber	bmRequestType, 714
CMcsUsbNet, 332	bRequest, 714
TxnTestMemoryReadAndCheck	wIndex, 714
CMcsUsbNet, 332	wLength, 714
TxnTestMemoryWrite	wValue, 715
CMcsUsbNet, 332	UsbTest
Туре	Mcs::Usb, 79
HeadStageIDType, 696	UsbVendorldEnumNet
TypeValue	Mcs::Usb, 90
HeadStageIDType, 696	User_ADC_0
	Mcs::Usb, 68
Unknown	User_ADC_1
HeadStageIDType, 694	Mcs::Usb, 68
Mcs::Usb, 52, 53, 65, 66, 87, 90	User_ADC_2
unknown	Mcs::Usb, 68
Mcs::Usb, 65	
UnknownDest	User_ADC_3
Mcs::Usb, 57	Mcs::Usb, 68
UnknownSpeed	User_ADC_4
Mcs::Usb, 70	Mcs::Usb, 68
Unlock	User_DAC_0
	Mcs::Usb, 69
Mcs::Usb, 76 UnlockPlateClamp	User_DAC_1
•	Mcs::Usb, 69
CMultiwellDeviceNet, 391	User_DAC_2
UnrolledAmplitude	Mcs::Usb, 69
CStimulusFunctionNet::StimulusDeviceDataAndUnro	USel DefinedName
713	HeadStageIDType, 696
UnrolledDuration	UssingChamber
CStimulusFunctionNet::StimulusDeviceDataAndUnro	
713	UssingClampModeEnumNet
UnrolledSync	Mcs::Usb, 90
CStimulusFunctionNet::StimulusDeviceDataAndUnro	ollessingaail
713	Mcs::Usb, 67
Unsigned_16bit	UssingUnitEnumNet
Mcs::Usb, 61	Mcs::Usb, 91
Unsigned_24bit	
Mcs::Usb, 61	Valid
Unsigned_32bit	HeadStageIDType, 696
Mcs::Usb, 61	ValidKey
UpdateChannelBlock	CMcsUsbNet, 332
CCMOSMeaDeviceNet, 115	Value
UpdateDisplay	CMcsUsbFactoryNet::TargetSelectionInfoNet, 714
CRoboDacqNet, 470	VendorldEnumNet
UpdateFirmware	Mcs::Usb, 91
CMcsUsbFactoryNet, 305–307	VendorInRequest
UpdateTransistorVoltages	CGenericDevelopDeviceNet, 167
CCMOSMea_FunctionNet, 112	VendorOutRequest
232323	. S S Sull loquoti

CGenericDevelopDeviceNet, 167	CWarnerUssingFunctionNet, 646
VirtualDevice_ContinousDacq	WaitTimer
CRoboDacqNet, 470	CRoboDeviceNet, 482
VirtualDevice_TableRun	Warner
CRoboDacqNet, 470	Mcs::Usb, 87
Volt	Warner_TEER_Machine
Mcs::Usb, 53, 91	Mcs::Usb, 79
	Warner_Ussing
BatteryState, 97	Mcs::Usb, 79
	WARNER_USSING_DEVICE
Mcs::Usb, 68	Mcs::Usb, 62
- -	Warner_Valve_Control
Mcs::Usb, 68	Mcs::Usb, 79
÷ ,	WARNER_VALVE_CONTROL_DEVICE
Mcs::Usb, 90	Mcs::Usb, 62
	WarnerUssingFunction
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	
	WClassicFunctionNet
VoltageRangeInMicroVolt	CMeaDeviceNet, 361
CMcsUsbDacqNet::CHWInfo::CVoltageRangeInfoNet	
614	Mcs::Usb, 78
	WholeCellPatch
VoltageResolutionInMicroVolt	Mcs::Usb, 82
	WholeCellPatchHeadstage
VoltageString	Mcs::Usb, 67
•	WidthFrequency
VOPSTimerSetResetTimes	CCreateFilterNet, 118
CCMOSMea_FunctionNet, 112	wIndex
W16lsW14	usbSetupPacket_t, 714
HeadStageIDType, 696	WirelessHeadStageAccDataRE1HS1
14/2/22	Mcs::Usb, 92
Mcs::Usb, 79	WirelessHeadStageAccDataRE1HS2
	Mcs::Usb, 92
Mcs::Usb, 91	WirelessHeadStageAccDataRE1HS3
MOLOO E N	Mcs::Usb, 92
CMeaDeviceNet, 361	WirelessHeadStageAccDataRE1HS4
14/04/00 Ott D 1 14/5	Mcs::Usb, 92
CurrentRangeInNanoAmp, 715	WirelessHeadStageAccDataRE2HS1
O 1D 11: 1A1 A 745	Mcs::Usb, 92
DACResolution, 715	WirelessHeadStageAccDataRE2HS2
Ti D 1 ii 1 N 0 1 7 7 5	Mcs::Usb, 92
VoltageRangeInMicroVolt, 715	WirelessHeadStageAccDataRE2HS3
V 15 B 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1	Mcs::Usb, 92
W2100DacqGroupChannelEnumNet	WirelessHeadStageAccDataRE2HS4
Manual July 04	Mcs::Usb, 92
W2100DigitalSourceEnumNet	WirelessHeadStageAnalogRE1HS1
Manual India 00	Mcs::Usb, 92
W2100IFB2	WirelessHeadStageAnalogRE1HS2
Manufile 07	Mcs::Usb, 92
W2100Interfaceboard	WirelessHeadStageAnalogRE1HS3
	Mcs::Usb, 92
W2100WirelessReceiver	WirelessHeadStageAnalogRE1HS4
M 111 07 00	Mcs::Usb, 92
W2100WirelessReceiverAnalog	WirelessHeadStageAnalogRE2HS1
M 111 07 00	Mcs::Usb, 92
WaitForAllChambers	WirelessHeadStageAnalogRE2HS2
OM-manullanian Franctica Net 040	Mcs::Usb, 92
WaitForChamber	WirelessHeadStageAnalogRE2HS3

Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageAnalogRE2HS4	WirelessHeadStageReservedBRE2HS1
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE1HS1	WirelessHeadStageReservedBRE2HS2
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE1HS2	WirelessHeadStageReservedBRE2HS3
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE1HS3	WirelessHeadStageReservedBRE2HS4
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE1HS4	WirelessHeadStageReservedCRE1HS1
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE2HS1	WirelessHeadStageReservedCRE1HS2
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE2HS2	WirelessHeadStageReservedCRE1HS3
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE2HS3	WirelessHeadStageReservedCRE1HS4
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageGyroDataRE2HS4	WirelessHeadStageReservedCRE2HS1
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageOptoStimCurrentRE1HS1	WirelessHeadStageReservedCRE2HS2
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageOptoStimCurrentRE1HS2	WirelessHeadStageReservedCRE2HS3
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageOptoStimCurrentRE1HS3	WirelessHeadStageReservedCRE2HS4
Mcs::Usb, 92	Mcs::Usb, 93
WirelessHeadStageOptoStimCurrentRE1HS4	WirelessHeadStageStatusRE1HS1
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageOptoStimCurrentRE2HS1	WirelessHeadStageStatusRE1HS2
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageOptoStimCurrentRE2HS2	WirelessHeadStageStatusRE1HS3
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageOptoStimCurrentRE2HS3	WirelessHeadStageStatusRE1HS4
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageOptoStimCurrentRE2HS4	WirelessHeadStageStatusRE2HS1
Mcs::Usb, 93	Mcs::Usb, 92
WirelessHeadStageReservedARE1HS1	WirelessHeadStageStatusRE2HS2
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageReservedARE1HS2	WirelessHeadStageStatusRE2HS3
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageReservedARE1HS3	WirelessHeadStageStatusRE2HS4
Mcs::Usb, 92	Mcs::Usb, 92
WirelessHeadStageReservedARE1HS4	WirelessTestAdapter
Mcs::Usb, 92	Mcs::Usb, 52
WirelessHeadStageReservedARE2HS1	wLength
Mcs::Usb, 92	usbSetupPacket_t, 714
WirelessHeadStageReservedARE2HS2	Work
Mcs::Usb, 92	Mcs::Usb, 94
WirelessHeadStageReservedARE2HS3	WPA16
Mcs::Usb, 93	Mcs::Usb, 79
WirelessHeadStageReservedARE2HS4	WPA32
Mcs::Usb, 93	Mcs::Usb, 79
WirelessHeadStageReservedBRE1HS1	WPA4
Mcs::Usb, 93	Mcs::Usb, 79
WirelessHeadStageReservedBRE1HS2	WPA8
Mcs::Usb, 93	Mcs::Usb, 79
WirelessHeadStageReservedBRE1HS3	WPAError_ScanningIsPending
Mcs::Usb, 93	CMcsUsbNet, 338
WirelessHeadStageReservedBRE1HS4	Write

```
CExternDTesterDeviceNet, 134
Write2
    CExternDTesterDeviceNet, 135
WriteEepromRegisterPreconfig
    CMcsUsbNet, 333
WritePipe
    CGenericDevelopDeviceNet, 167
WriteRegister
    CMcsUsbNet, 333
WriteRegister32
    CMcsUsbNet, 333
WriteRegisterArray
    CMcsUsbNet, 333
WriteRegisterTimeSlot
    CMcsUsbNet, 334
WriteRegisterValue
    CMcsUsbNet, 334
WriteUARTData
    CLIH3DeviceNet, 204
wValue
    usbSetupPacket_t, 715
WvcDisplayModeEnumNet
    Mcs::Usb, 94
WvcValveModeEnumNet
    Mcs::Usb, 94
Zero
    Mcs::Usb, 57, 71, 73, 74, 83, 86, 93
```